
UV and IR Spectra

Pharmaceutical Substances (UV and IR) and Pharmaceutical and Cosmetic Excipients (IR)

Introduction in English and German

H.-W. Dibbern R. M. Müller E. Wirbitzki (Eds.)



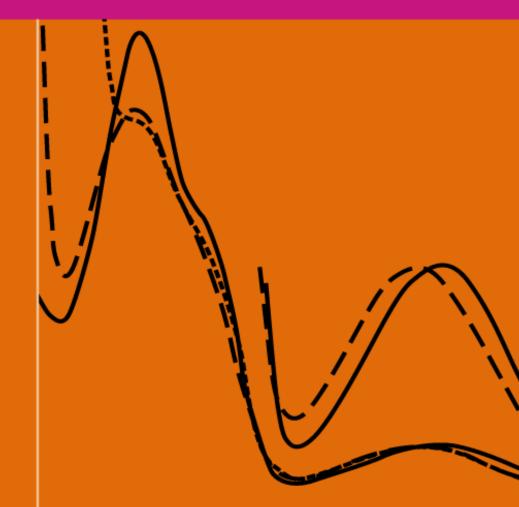

ISBN 3-87193-238-8 © 2002 ECV · EDITIO CANTOR VERLAG

General

- Functions of the CD
- Funktionalität der CD
- Preface
- Vorwort
- Introduction (Pharmaceutical substances)
- Introduction (Excipients)
- Einleitung (Wirkstoffe)
- Einleitung (Hilfsstoffe)
- Further information
- Weitere Informationen

Funktionalität der CD

Diese CD-ROM bietet Ihnen diverse Arbeitsmöglichkeiten mit insgesamt 1673 Spektren.

Die allgemeinen Funktionalitäten sind identisch mit denen des Acrobat Reader. Eine ausführliche Dokumentation dafür finden Sie unter dem Menü Hilfe im Punkt Reader Hilfe.

In den **Allgemeinen Informationen** sind genaue Verlagsangaben, Lizenzbedingungen und Hinweise für weitere Informationen zum Verlagsprogramm enthalten.

Die Texte zum **Vorwort**, zu den Allgemeinen Informationen, zu den **Benutzerhinweisen** und zur **Einleitung** sind in englischer und in deutscher Sprache aufrufbar. Auch die **Spektrokognosie** ist zweisprachig aufgenommen. Alle weiteren Dokumente der Spektrensammlung sind in englischer Sprache abgefaßt.

Die **Spektrensammlung** selbst ist aufgeteilt in die **Hauptgruppen** *Wirkstoffe* (dazu auch: Indikationen) **und** *Hilfsstoffe*. Sie können Texte bzw. auch einzelne Graphikelemente markieren und kopieren für die Weiterverwendung in anderen Programmen. Ebenso können Sie ganze Blätter ausdrucken.

In der **Specfinder**-Tabelle können Sie nach UV-Absorptionsmaxima suchen. Geben Sie dazu in der Suchfunktion des Acrobat Reader die gewünschte Zahl + Angabe nm ein (ohne Leerschritt; z.B.: 348nm). Das gefundene Spektrum läßt sich durch Anklicken der entsprechenden Spektren-Nummer aufrufen. (Für ausführlichere Angaben vgl. die Einleitung.)

Die **weiteren Tabellen** enthalten Dokumente für "Wirkstoffe, die nur in alkalisch-(wäßriger) Lösung ein differenziertes UV-Spektrum zeigen oder in Methanol nicht ausreichend löslich sind", für "Wirkstoffe ohne Absorptionsmaximum" und für "Wirkstoffe ohne analytisch verwertbare UV-Absorption > 220 nm". Hier können Sie nur recherchieren, eine Verlinkung ist nicht vorhanden.

Die Bildschirmdarstellung der Spektrenblätter sowie der weiteren Dokumente richtet sich zunächst nach den Komponenten Ihres PC sowie nach den Einstellungen des Acrobat Reader. Vergrößerungen, ggf. auch Verkleinerungen können Sie mit den Funktionalitäten dieses Programmes selbst bestimmen.

Für Fragen, Kritik oder Anregungen zu diesem Produkt nutzen Sie bitte die folgenden Kontaktmöglichkeiten:

e-mail: support@ecv.de

oder

Fax +49 (0) 75 25-94 01 80

Functions of the CD

This CD-ROM with its 1,673 spectra provides various working options.

The general functions are identical with those of Acrobat Reader. For detailed information please refer to *Reader Help* under the menu *Help*.

Under **Further information** you will find details about the publication and the publishing company, license conditions and indications concerning the publisher's products.

Preface, Further information, Functions of the CD, and **Introductions** are available both in English and German, and so is the **Spectrocognosia.** All other documents of the collection of spectra are in English.

The **collection of spectra** is divided into **two major groups: pharmaceutical substances** (including indications) and **excipients.** You can highlight texts or single graphs and copy them for further use in other applications. You can also print whole sheets.

In the **Specfinder table** you can search for UV absorption maximum values. This is possible by entering the particular figure + nm (without blank; e. g. 348nm) in the search function of Acrobat Reader. The found spectrum can be called up by clicking on the respective spectrum number. (See the Introduction for more detailed information.)

The section **Further tables** contains documents for "pharmaceutical substances showing a well distinct UV spectrum in alkaline solution or being unsoluble in methanol", "pharmaceutical substances showing no distinct maximum of absorption", and "pharmaceutical substances without useful UV absorption > 220 nm". There you can only search, a link is not available.

The display of the spectra sheets or other documents on the screen depends both on the components of your PC and the settings of Acrobat Reader. Magnification resp. diminution is possible using the functions of this application.

In case of any questions, criticism or suggestions, please contact:

e-mail: support@ecv.de

or

Fax +49 (0)7525-940 180

Vorwort

Unter dem Titel "UV- und IR-Spektren wichtiger pharmazeutischer Wirkstoffe" ist 1978 erstmals eine Loseblatt-Sammlung von Spektren erschienen, die durch mehrere Nachlieferungen bis 1998 erweitert wurde. 1987 wurde eine zweite Sammlung mit "IR-Spektren von pharmazeutischen und kosmetischen Hilfsstoffen" herausgegeben.

Überlegungen mit Blick auf eine Fortführung der beiden Werke, die mittlerweile einen Umfang von fünf DIN-A 4-Ordnern angenommen hatten, mußten zu einer neuen Konzeption führen, die dennoch inhaltlich den bisherigen Sammmlungen Rechnung trägt. Das Ergebnis konnte nur eine elektronische Version darstellen, die beide Werke aufnimmt und bei der die Texte in englischer Sprache abgefaßt sind (Hinweise für die Nutzung z. B. sind indessen zweisprachig Englisch/Deutsch).

Die Spektrensammlungen liegen nun in Form dieser CD-ROM-Ausgabe vor, bei der Spektren von 55 weiteren Wirksubstanzen Aufnahme gefunden haben.

Die Absicht der Herausgeber und des Verlages ist, mit den unzweifelhaften Vorteilen eines elektronischen Mediums und dem damit verbundenen zusätzlichen Komfort den Benutzern den Gebrauch der Sammlung erheblich zu erleichtern. Dabei bleibt allerdings wesentliche Voraussetzung für die Sinnhaftigkeit einer solchen Sammlung, was bereits in der Einleitung zur früheren Ausgabe der UV- und IR-Spektren ausgeführt wurde – daß nämlich für eine möglichst sichere Zuordnung von Spektren zu bestimmten pharmazeutisch wichtigen Wirkstoffen neben einigen exakten Daten ganz besonders auch die "individuelle" Form der Spektren – gleichsam einem 'Fahndungsphoto' – von entscheidender Bedeutung für die Identifizierung oder doch engmaschige Zuordnung ist.

In dem Verzeichnis und bei den Spektren sind die gebräuchlichen Kurzbezeichnungen INN (International Nonproprietary Names), INNv (vorgeschlagene internationale Kurzbezeichnungen) und sonstige übliche Bezeichnungen aufgeführt. Zusätzlich sind bei den Spektren aller Wirkstoffe (Ausnahme: komplexe Moleküle, Polymere) die chemischen Strukturformeln und die relativen Molekülmassen angegeben. Auf die Bezeichnung über die wissenschaftliche Nomenklatur wurde daher verzichtet.

Das System der Numerierung nach Indikationsgruppen (bei den Wirkstoffen) bzw. nach Stoffklassen (bei den Hilfsstoffen) wurde belassen. Zusätzlich sind alle Wirk- und Hilfsstoffe in einem alphabetischen Verzeichnis aufgelistet. (Vgl. hierzu die Hinweise in der Einleitung.)

Das UV-Specfinder-System wurde als bewährtes Hilfsmittel in die Sammlung integriert und ist jetzt mit den Vorteilen einer elektronischen Version nutzbar.

Im Verlaufe der sich über viele Jahre erstreckenden ständigen Erweiterungen der Sammlungen erfolgten zwangsläufig auch immer wieder Detailänderungen auf den Spektrenblättern. Somit kann sich beim Zusammenfügen aller Spektren zwangsläufig kein absolut einheitliches Bild in der optischen Präsentation ergeben. Der Nutzwert indessen ist in keiner Weise beeinträchtigt.

Allen Personen, die in Unternehmen der Pharma-Branche (insbesondere in den Bereichen der Analytik, Forschung und Entwicklung, Qualitätskontrolle u. ä.), in der Zulieferindustrie, in Untersuchungsämtern, Überwachungsbehörden, in Universitäten, in der Gerichtsmedizin etc. tätig sind, wünschen wir einen nutzbringenden und hilfreichen Einsatz unserer Spektrensammlung in der neuen, komfortablen Form der CD-ROM-Version.

Frankfurt/Main und Leofels, September 2002

H.-W. Dibbern, R.M. Müller und E. Wirbitzki

Preface

In 1978, a loose-leaf spectra collection entitled "UV and IR Spectra of Some Important Drugs" first appeared, which until 1998 was expanded by several supplements. In 1987, a second collection was published: "IR Spectra of Pharmaceutical and Cosmetic Auxiliary Substances".

The idea of continuing both works which meanwhile comprise five files necessarily led to a new conception, however, taking into account the contents of the previous collections. The result could only be an electronic version including either work and with the texts being written in English (e. g. user information being bilingual English/German).

The spectra collections are available on this CD-ROM now including the spectra of 55 additional pharmaceutical substances.

It was the intention of both the editors and the publisher to enable the user to more easily use the collection with the aid of the indubitable advantages of an electronic medium and the conveniences associated therewith. However, as already emphasized in the introduction of the previous edition of UV and IR Spectra there is one essential prerequisite for the meaning of such a collection: Apart from some exact data it is in particular the "individual" form of a spectrum – the "wanted" photo so to speak – which is of utmost importance for the identification or at least most probable assignment in order to be able to reliably assign a spectrum to a particular pharmaceutical substance.

The usual names such as INN (International Nonproprietary Names), INNv (proposed International Nonproprietary Names) and other customary names are given in the indexes. In addition, with the spectra of all pharmaceutical substances (with the exception of complex molecules and polymers) the chemical structural formulas and the relative molecular weights are indicated. Therefore the designation by scientific nomenclature has been dispensed with.

The numbering system by therapeutic areas (for the pharmaceutical substances) and by classes of substances is the same as in the previous edition. In addition all pharmaceutical substances and excipients are listed in an alphabetical index. (See the respective information in the Introduction.)

The proven UV Specfinder system has been integrated in the collection offering now the advantages of an electronic version.

In the course of the continuous expansions of the collections over many years some changes in the spectra sheets have been made repeatedly. Thus putting together all spectra the appearance inevitably may not be uniform. However, that does not impair the utility value.

We do hope that all those active in the pharmaceutical industry (in particular in analytics, research and development, quality control or related areas), ancillary industry, public health authorities or public health labs, regulatory authorities, universities, forensic medicine etc. will benefit from the use of this new and convenient collection of spectra on CD-ROM.

Frankfurt/Main and Leofels, September 2002

H.-W. Dibbern, R.M. Müller and E. Wirbitzki

Einleitung (Hilfsstoffe)

Außer den Wirkstoffen haben auch die für Fertigarzneimittel verwendeten Hilfsstoffe eine – zuweilen recht wesentliche – Bedeutung für den therapeutischen Wert eines Arzneimittels. U. a. können sie für die biologische Verfügbarkeit und für die Haltbarkeit entscheidend mitbestimmend sein.

Da UV-Spektren in dieser Stoffklasse in aller Regel unergiebig sind, haben wir uns in der Sammlung der Hilfsstoffe auf IR-Spektren beschränkt. Diejenigen Hilfsstoffe, die auswertbare UV-Spektren erbringen, sind in der Sammlung »UV- und IR-Spektren wichtiger pharmazeutischer Wirkstoffe« bereits aufgeführt. Tritt bei der Analyse eine nicht als Wirkstoff deklarierte Substanz auf, die eine auswertbare UV-Absorption zeigt, dann sollte unsere Sammlung Wirkstoffspektren herangezogen werden. Selbstverständlich ist auch die Abwesenheit einer UV-Absorption ein analytisch brauchbarer Hinweis. Für die weitere Kennzeichnung eignet sich aber ihr IR-Spektrum, allerdings mit folgenden Einschränkungen:

Pharmazeutische und kosmetische Hilfsstoffe sind zwar vorwiegend chemisch einigermaßen definierte Verbindungen, doch gilt dies nur in Annäherung, da viele von ihnen Naturstoffe sind oder von solchen abstammen. Daher muß mit Abweichungen von der strengen Definition – z. B. bei der Kettenlänge – häufig gerechnet werden. Dies hat aber in der Regel praktisch keinen Einfluß auf das IR-Spektrum.

Es war unser Bestreben, jeweils für bestimmte Stoffklassen typische und häufig gebrauchte Vertreter auszuwählen. Diese können dann als Repräsentanten für ähnliche Verbindungen herangezogen werden. Dies ist auch der Grund, weshalb die Zahl der hier aufgenommenen Stoffe im Vergleich zu Listen über Hilfsstoffe, z. B. »Fiedler · Lexikon der Hilfsstoffe« oder »Index der Hilfsstoffe« von H. P. Fiedler, vergleichsweise klein ist. Es wäre nicht sinnvoll gewesen, zahlreiche sehr ähnliche Substanzen – die alle praktisch identische IR-Spektren zeigen – aufzunehmen.

Aus unserer Kenntnis über sehr zahlreiche Hilfsstoffspektren haben wir folgende Einteilung in Stoffklassen aufgestellt:

- 1. Organische Säuren und deren Salze; Aminosäuren
- 2. Organische Basen und deren Salze
- 3. Kürzer- und mittelkettige Ester
- 4. Fettsäuren und deren Salze, Fettsäurederivate, Fette
- 5. Fettalkohole, Fettalkoholderivate
- 6. Kohlenwasserstoffe
- 7. Wachse und deren Derivate
- 8. Kohlehydrate
- 8.1. Leicht lösliche Zucker und Zuckeralkohole
- 8.2. Wenig lösliche oder nur quellbare Cellulose- und Stärkederivate; andere Polysaccharide
- 9. Anorganische Hilfsstoffe, wasserunlöslich
- 10. Silicone
- 11. Synthetische Polymere
- 12. Sorbitanderivate
- 13. Glykole und andere mehrwertige Alkohole und deren Derivate (außer Fette)
- 14. Sonstige Naturstoffe bzw. Naturstoffgemische

Innerhalb dieser Stoffklassen wurden diejenigen Einzelstoffe aufgenommen, die im IR-Spektrum noch eine deutliche Differenzierung ermöglichen. Wenn in der Sammlung aber doch zwei Substanzen mit sehr ähnlichen IR-Spektren auftreten, so liegt der Grund darin, daß eine derartige Ähnlichkeit oft nicht zu erwarten war.

Einleitung (Hilfsstoffe) Fortsetzung

Überhaupt muß an dieser Stelle festgestellt werden, daß – anders als bei Wirkstoffen – eine absolute Identifizierung über das IR-Spektrum nur in wenigen Fällen erreichbar ist. Für die Praxis erscheint es uns aber wichtig, zunächst wenigstens eine Zuordnung zu einer der o. g. Stoffklassen zu erreichen. Häufig ist dann auch eine nähere Zuordnung zu einem der in der Sammlung aufgeführten Stoffe möglich. Erleichtert wird dies durch unsere Einteilung in die Stoffklassen, bei denen zum Teil auch Löslichkeitseigenschaften mit berücksichtigt worden sind: Ist eine fragliche Substanz einer solchen Stoffklasse zugeordnet worden, so ist es verhältnismäßig leicht, einer Vorklassifizierung nahe zu kommen, indem die Spektren der in dieser Stoffgruppe aufgeführten Substanzen mit dem fraglichen Spektrum verglichen werden. Eine eindeutige Identifizierung eines auf diese Weise vorklassifizierten Stoffes dürfte anschließend mit anderen Mitteln (GC, MS, NMR) relativ leicht möglich sein.

Wegen der größtenteils nur wenig differenzierten Spektren erschien uns der Aufbau eines Specfinder-Systems hier nicht aussichtsreich. In Anbetracht der kleinen Zahl beispielhafter Spektren ist es nach unserer Erfahrung aber relativ leicht möglich, mittels Durchsehen der Sammlung zu einer Zuordnung zu kommen.

In den meisten Fällen haben wir für die Zuordnung und für die Auffindung einer Substanz die englische chemische Bezeichnung verwendet; nicht möglich ist dies naturgemäß bei »Hilfsstoffen«, die aus mehreren Komponenten aus zum Teil sehr verschiedenen chemischen Stoffklassen zusammengesetzt sind. Hier blieb uns keine andere Wahl, als allein die Markenbezeichnungen zu benützen. Die Nennung eines solchen Warenzeichens besagt nicht, daß nicht auch noch andere Präparate gleicher oder ähnlicher Zusammensetzung auf dem Markt sind.

Die Nummern gemäß Stoffklasseneinteilung befinden sich auf dem Spektrum in dessen rechten unteren Ecke.

Spektrenaufnahme: Die IR-Spektren wurden mit einem Perkin-Elmer Spektrophotometer Modell 298 mit einer Geschwindigkeit von 12 Minuten pro Spektrum aufgenommen. Spektrale Auflösung des Photometers: 3 cm⁻¹ bei 1100 cm⁻¹.

Angaben zur Probenpräparation befinden sich auf jedem Spektrenblatt in der rechten oberen Ecke.

Die Darstellung der IR-Spektren erfolgt wie allgemein üblich: Ordinate in Prozent Transmission, Abszisse als Wellenzahl bzw. Wellenlänge.

Einleitung (Wirkstoffe)

Zweck und Gebrauch der Spektrensammlung

Die Chance, bei einer Fahndung einen Unbekannten nur auf der Basis einiger Einzelangaben zu erkennen und aufzufinden, ist bekanntlich recht gering: Mitteilungen über Größe, Gewicht, Alter, Farbe der Augen und Haartracht bringen nur eine bescheidene Einengung aus einer großen Zahl von Individuen. Erst die Summe zahlreicher Merkmale – die großenteils unbewußt und kaum definierbar wahrgenommen werden – macht ein Erkennen möglich. Aus diesem Grund stellt eine anschauliche Abbildung – womöglich zusammen mit exakten Daten – eine ungleich bessere Information dar.

Ähnlich ist es auch in der Spektroskopie: UV- und IR-Spektren einer unbekannten Substanz auf »theoretischer« Grundlage zu bewerten, ist schwierig und erlaubt auch dem Spezialisten nur selten eine sichere Zuordnung.

Spektrophotometrische Untersuchungen gründen ihre Auswertung – qualitativer wie quantitativer Zielsetzung – auf Vergleiche.

Absorptionsmaxima und -minima oder spezifische Extinktionen im UV lassen sich in Zahlenwerten ausdrücken. Diese Daten sind wertvolle Wegweiser zur Auffindung und werden auch in den folgenden Tabellen verwendet; aber auch die Breite einer Absorptionsbande oder eines Minimumbereichs, die Steilheit, die Wendepunkte oder die geringfügigen Schultern einer Kurve sind für die Zuordnung oft unverzichtbare Merkmale eines UV-Spektrums – wie eben die Form der Nase, der Schnitt des Mundes, Rundungen und Falten eigentlich erst die Physiognomie ausmachen!

Aus dieser Anschauung heraus haben wir im Verlauf der letzten Jahre eine UV-Spektrensammlung wichtiger Arzneistoffe aufgebaut, die uns stets einen schnellen Zugriff zu wichtigem Vergleichsmaterial ermöglicht und damit die analytische (Fahndungs-)Arbeit wesentlich erleichtert. Später wurde die Sammlung durch Hinzufügung der entsprechenden IR-Spektren ergänzt.

Manche Kollegen aus anderen Arbeitsgruppen haben in der Vergangenheit unsere Sammlung gerne mitbenutzt. So lag es nahe, eine geeignete Auswahl dieser Sammlung einem größeren Kreis von Benutzern zugänglich zu machen.

Was bedeutet nun »geeignete Auswahl«? Es gibt einige tausend Wirkstoffe, die in pharmazeutischen Stofflisten (aller Art) aufgeführt werden; und es gibt eine kleinere Anzahl von Substanzen, die einen breiteren Eingang in der Arzneimitteltherapie gefunden haben. Eine gewisse Richtschnur hierfür sind die modernen Pharmakopöen; nach diesen haben wir uns im großen und ganzen orientiert. Allerdings haben wir einige Arzneibuchpräparate weggelassen, dafür aber eine größere Anzahl von neueren Wirkstoffen aufgenommen, die in Pharmakopöen nicht enthalten sind, die uns aber aus verschiedenen Gründen wichtig oder interessant erschienen. Daß eine derartige Auswahl nicht die uneingeschränkte Zustimmung aller Benutzer finden wird, ist wohl unvermeidlich. Es ist auch selbstverständlich, daß eine derartige Auswahl immer wieder der Ergänzung und Aktualisierung bedarf. Anregungen für eine Erweiterung der Sammlung nehmen wir dankbar auf.

Hinsichtlich der Nomenklatur sind wir für die Spektren von den nicht wortgeschützten Kurzbezeichnungen ausgegangen, und zwar, soweit vorhanden, von den empfohlenen Freinamen (INN). Diese Handhabung erschien uns konsequent, weil die Spektren sich auf die Wirksubstanzen beziehen, zu deren Kennzeichnung die Kurzbezeichnungen dienen, während die Zubereitungen in Gestalt der pharmazeutischen Spezialitäten mit Warenzeichen bezeichnet werden.

Für die Darstellung der UV-Spektren wurde als Ordinate die Absorption E_{1cm} der jeweils gemessenen Lösung (mit Angabe der Konzentration) gewählt. Dies scheint uns für den praktischen Gebrauch des Analytikers, vor allem für die quantitative Auswertung, am günstigsten. Zweifellos hat auch die sonst vielfach übliche Darstellung der Ordinate als log ϵ gewisse Vorzüge, doch haben wir die Wiedergabe der Spektren, »wie sie anfallen«, aus Gründen der praktischen Anwendung bewußt bevorzugt.

Einleitung (Wirkstoffe) Fortsetzung

Die Darstellung der IR-Spektren erfolgt, wie heute allgemein üblich: Ordinate in % Transmission, Abszisse als Wellenzahl bzw. Wellenlänge. Wenn nicht anders vermerkt, wurden hier KBr-Preßlinge 1,5 mg + 300 mg eingesetzt.

Nach unseren Erfahrungen kann die Spektrensammlung zur Lösung verschiedener Aufgaben behilflich sein:

Sie kann zur Identifizierung oder Zuordnung unbekannter Wirkstoffe beitragen. Selbstverständlich wird man heute in größeren Laboratorien derartige Probleme unter Mitheranziehung von NMR- und Massenspektroskopie lösen. Immerhin gelingt eine Zuordnung – wenigstens zu einer Substanzklasse – in vielen Fällen auch schon aus UV-Spektren – sofern Vergleichsmaterial bzw. Vergleichsspektren zur Verfügung stehen. Häufig wird dabei ein bestimmtes Indikationsgebiet als wahrscheinlich vorausgesetzt werden können, wie z. B. »Analgetikum«, »Lokalanästhetikum«, »Psychopharmakon« oder ähnliches. In solchen Fällen kann unter »Indications« der primär interessierende Bereich ausgewählt werden. Anschließend können die in Frage kommenden Spektrenblätter aufgerufen werden. Gegebenenfalls ist der betreffende Wirkstoff auch in zwei oder drei Indikationsgruppen zu finden.

Für eine zuverlässige Identifizierung reichen UV-Absorptionsdaten allein meist nicht aus. Deshalb müssen auch IR-Spektren mit aufgenommen werden.

Für quantitative Auswertungen sind den UV-Spektren die Absorptionsdaten bei Verwendung verschiedener Lösungsmittel in Form der E 1%/1cm und der ϵ -Werte beigefügt. Die in der Tabelle jeweils angegebenen Werte wurden von im Handel verfügbaren Reinsubstanzen erhalten und auf mehreren Photometern überprüft. Daß bei diesen Werten kleine Differenzen mit den Befunden anderer Untersucher möglich sind, wird jedem Kenner der Materie klar sein.

Die Wiedergabe der Strukturformel soll dem Benutzer dieser Sammlung die Zusammenhänge zwischen chemischer Konstitution und Absorptionseigenschaften nahebringen mit dem Ziel, für sich selbst eine Spektrokognosie zu entwickeln, die ihn befähigt, ein bekanntes Spektrenbild zu erkennen und die Spektren unbekannter Komponenten – wenigstens mit einer guten Wahrscheinlichkeit – einem Verbindungstyp zuzuordnen oder gewisse funktionelle Strukturelemente zu erkennen – so etwa phenolische OH- oder aromatische Aminogruppen, Pyrazolone, Xanthine, Phenothiazine usw. Wir haben einige Beispiele zu einer solchen Arzneimittel-UV-Spektrokognosie zusammengestellt und dabei gut erkennbare Zusammenhänge zwischen Struktur und Spektrum aufgezeigt. Diese Beispiele sollen als Wegweiser für die Praxis dienen; sie sind aber auch als Anregung gedacht, einen derartigen »Erkennungsdienst« für spezielle Stoffklassen selbst weiterzuführen.


Ganz bewußt stellen wir für die Aufgabe der Zuordnung das übersichtlichere UV-Spektrum an die erste Stelle.

Selbstverständlich sind solche Charakterisierungen und Zuordnungen nicht überall gleich gut möglich, und häufig reicht die Information UV-Spektrum allein für eine Kennzeichnung nicht aus. Immerhin sollte es aber doch in einigen Fällen möglich sein, mit Hilfe solcher spektrokognostischen Abbildungen und Hinweise manche unbekannte Substanz, die in der Sammlung nicht aufgeführt ist, schon aufgrund ihres UV-Spektrums relativ schnell zu einer durch ihre Chromophore geprägten chemischen (und damit öfters auch einer therapeutischen) Gruppe zuzuordnen. So liegt es z. B. nahe, bei einem Pyrazolon-Derivat auf die Gruppe der Analgetika/Antipyretika zu schließen. Umgekehrt wird man bei Verdacht auf ein Lokalanästhetikum an p-Aminobenzoesäureester oder an basisch substituierte Anilide denken und das Spektrum unter diesem Aspekt beurteilen.

Liegt keinerlei Hinweis auf eine chemische Gruppe und auch keine Indikationszuordnung vor, so kann man den Weg über die Liste der Absorptionsdaten gehen, um die Zuordnung eines Spektrums zu versuchen. Eine Vollidentifizierung unbekannter Stoffe aus UV-Absorptionsdaten allein ist allerdings selten möglich. Innerhalb der hier vorliegenden Sammlung und unter Heranziehung der Spektren, auf die die Specfinder-Liste hinführt, gelingt aber eben doch oft eine Zuordnung oder Kennzeichnung. Zur eigentlichen Identifizierung einer Substanz ist dann das IR-Spektrum mit heranzuziehen.

Als erstes Zuordnungsprinzip (I) benutzen wir in unserem UV-Specfinder jeweils das (deutlich ausgeprägte) längstwellige UV-Absorptionsmaximum in methanolischer Lösung. Dieses kann auch das einzige oder das stärkste Maximum des gesamten hier aufgenommenen Bereichs sein. Das längerwellige Maximum kann aber auch – gegenüber kürzerwelligen – eine geringere spezifische Absorption besitzen; dennoch soll das Maximum im Längerwelligen als Zuordnungsmerkmal I dienen. Sind in der Liste zwei (oder drei) Maxima aufgeführt und ist eines der unter II genannten Maxima stärker als das langwellige »Leit-Maximum« I, so ist das stärkste Maximum durch Unterstreichung gekennzeichnet. Die Frage, ob das längstwellige Maximum »hinreichend deutlich« ist, mag in einzelnen Fällen strittig sein. In diesen Fällen suche man stärkere Maxima bei kürzerer Wellenlänge.

Liegen mehrere Feinbanden in einer »Gruppe« relativ eng zusammen, so gilt hier die stärkste Bande dieses Absorptionsbereichs als maßgebend für die Zuordnung. Beispiel: Tolbutamide Sp.-Nr. 508. Erstes Zuordnungsprinzip I ist die stärkste Bande des langwelligen Bereichs 263 nm (nicht 273), zweites Zuordnungsprinzip II ist die – wesentlich stärkere – kurzwellige Bande in Methanol bei <u>228</u> nm.

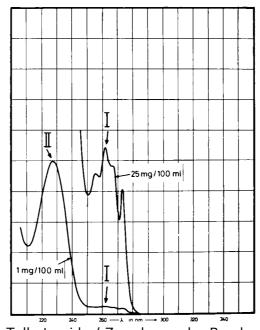

Tolbutamide / Zuordnung der Banden

Hierzu ein weiteres Bestimmungsbeispiel nach der Tabelle: Langwelliges ausgeprägtes Leit-Maximum I bei 330 nm, Maxima bei 254 und 219 nm = Chloroquine phosphate (Sp.-Nr. 908).

Bei einigen Wirkstoffen (z. B. sehr deutlich bei solchen mit phenolischen OH- oder aromatischen Aminogruppen) kommt es – im Vergleich zu Methanol – im wäßrig alkalischen oder sauren Milieu infolge Änderung des Dissoziationsgrades zu mehr oder weniger starken Deformationen des Spektrums. Solche Veränderungen sind in der Spalte IV der Specfinder-Tabelle durch folgende Symbole gekennzeichnet:

- † Erhöhung des Absorptionsmaximums
- ↓ Erniedrigung des Absorptionsmaximums
- 1 Erhöhung eines Maximums, Erniedrigung eines anderen
- → bathochrome Verschiebung
- ← hypsochrome Verschiebung
- 7 Erhöhung und bathochrome Verschiebung
- Erniedrigung und hypsochrome Verschiebung
- X weitgehende Veränderung des Spektrums
- X 283 weitgehende Veränderung des Spektrums und Auftreten eines neuen Maximums bei 283 nm
- →Z zunehmende Zersetzung

Einleitung (Wirkstoffe) Fortsetzung

Bei einigen Substanzen erhält man nur in verdünntem Alkali (0,1 M Natriumhydroxid) ausgeprägte Spektren. In solchen Fällen ist auch nur dieses Spektrum angegeben. Für diese Spektren wurde eine eigene Tabelle aufgenommen. Andere Wirkstoffe sind in wässrigen Lösungsmitteln zu wenig löslich, hier liegt nur ein Spektrum in Methanol vor.

Da die Angaben über UV-Absorptionsmaxima auf 1 nm genau naturgemäß mit einer gewissen Unsicherheit behaftet sind, sollten in jedem Fall auch die Werte bis \pm 2 nm der Liste bei der Zuordnung mit in Betracht gezogen werden.

Fast in allen Fällen werden auch quantitative Absorptionsdaten (E 1%/1cm) für die Zuordnung eines UV-Spektrums mit herangezogen. Daher muß hier von bekannten Konzentrationen ausgegangen werden.

Bis auf einige Ausnahmen sehr ähnlicher Verbindungen gelingt mit diesem zusammengesetzten Schlüssel eine durchaus brauchbare Zuordnung unter den hier aufgeführten Substanzen. Für die eigentliche Identifizierung stehen dann die IR-Spektren der in Betracht zu ziehenden Verbindungen zur Verfügung.

Ein weiterer wichtiger Anwendungsbereich der Sammlung liegt im Nachschlagen von Absorptionsdaten bekannter Arzneistoffe als Basisinformation für analytische Aufgaben, z. B. zur Beantwortung folgender Fragen:

Ist die Bestimmung des Stoffes X über die UV-Absorption überhaupt möglich, bzw. unter den gegebenen Umständen aussichtsreich?

oder:

Können die Stoffe X und Y (in den vorliegenden Konzentrationsverhältnissen) direkt oder über eine Simultanmessung bei 2 oder 3 Wellenlängen oder bei verschiedenen pH-Werten spektrophotometrisch bestimmt werden, oder ist eine vorhergehende (z. B. chromatographische) Trennung erforderlich oder doch wünschenswert?

Dem erfahrenen – aber auch dem angehenden – Analytiker werden diese Hinweise genügen, um auch noch weitere Anwendungsmöglichkeiten für die vorliegende Sammlung zu erkennen.

Wir wünschen allen Benutzern einen erfolgreichen Gebrauch. Für Anregungen aus der Praxis bezüglich notwendiger Änderungen oder wünschenswerter Ergänzungen sind wir stets dankbar.

Für technische Mitarbeit danken wir den Herren R. Merkel, M. Petry und W. Weiss.

Introduction (Excipients)

In addition to active pharmaceutical ingredients excipients used for the manufacture of finished drugs have – sometimes essential – influence on the therapeutic value of a drug. They may, e. g., be a contributory determinant of bioavailability and shelf-life.

As UV spectra for this class of substances are poor as a rule, in the collection we have confined ourselves to the IR spectra. Those excipients yielding evaluable UV spectra are already included in the collection of »UV and IR Spectra of Some Important Drugs«. If during analysis a compound not indicated as an active ingredient shows up, which presents evaluable UV absorption, our collection of drug spectra should be consulted. Of course also the absence of UV absorption is an analytically useful hint. For further characterization the IR spectrum is appropriate, however, with the following restrictions:

Pharmaceutical and cosmetic excipients are mostly compounds, which chemically are defined to a certain degree, but this holds true in approximation only, as many of them are naturally occurring substances or derived from such compounds. Therefore often deviations from the strict definition – e. g. regarding the chain length – have to be taken into account. However, this virtually is without impact on the IR spectrum.

We endeavoured to select typical and often used representative substances for particular classes of compounds. These then can be used as representatives for similar compounds. This also is the reason why the number of substances covered here is rather low compared with some lists of excipients such as, e. g., »Fiedler · Encyclopedia of Excipients« or »Index der Hilfsstoffe/Index of Auxiliary Substances« by H. P. Fiedler. It would not have been meaningful to include numerous, very similar substances – all showing virtually identical IR spectra.

According to our knowledge of a great number of spectra of excipients we arranged the division into classes of compounds as follows:

- 1. Organic acids and their salts; amino acids
- 2. Organic bases and their salts
- 3. Short and medium chain length esters
- 4. Fatty acids and their salts, fatty acid derivatives, fats
- 5. Fatty alcohols, fatty alcohol derivatives
- 6. Hydrocarbons
- 7. Waxes and their derivatives
- 8. Carbohydrates
- 8.1. Readily soluble sugars and sugar alcohols
- 8.2. Sparingly soluble or only swellable cellulose and starch derivatives; other polysaccharides
- 9. Inorganic excipients, water insoluble
- 10. Silicones
- 11. Synthetic polymers
- 12. Sorbitan derivatives
- 13. Glycols and other polyvalent alcohols and their derivatives (except fats)
- 14. Other naturally occurring compounds and mixtures of naturally occurring compounds

Within these classes of substances those single compounds have been included, which enable a distinct differentiation in the IR spectrum. However, if within this collection two substances with very similar IR spectra occur, the reason for this is that such a similarity often could not be expected.

In this connection it has to be stated anyway that – unlike for drugs – an absolute identification via the IR spectrum is achieved in few cases only. For the practical work, however, it seems important to us to achieve an assignment to one of the above mentioned classes of substances. In many cases a more precise assignment to a compound included in the collection is thus possible. This is facilitated by our division into classes of substances, partly also considering the solubility properties: Has a compound in question been assigned to such a class, it is comparatively easy to approach a precharacterization by comparing the spectra of the compounds included in this class of substances with the spectrum in question. A clear-cut identification of a thus pre-characterized compound surely is easily possible then using other methods (GC, MS, NMR).

Introduction (Excipients) continued

Due to the mostly poorly differentiated spectra putting up a specfinder system did not seem to be promising here. In view of the small number of exemplary spectra, in our experience it is easily possible to get to an assignment by searching the collection.

In most cases, for the assignment and discovery of a substance we have indicated the chemical name in English; naturally this is not possible with the »excipients«, which are composed of several components from partly very different chemical classes. Here we had no other choice than to use the trademarks. Indication of such a trademark does not imply that no other preparations of identical or similar composition are on the market.

The numbers according to the division into classes of substances are given on the spectrum on the lower right.

Taking the spectra: The IR spectra have been taken with a Perkin-Elmer Spectrophotometer Model 298 at a speed of 12 min per spectrum. Spectral resolution of the photometer was 3 cm⁻¹ with 1100 cm⁻¹.

Details for the sample preparation are given on each spectrum sheet on the upper right.

Presentation of the IR spectra is in accordance with common practice: ordinate in percent transmission, abscissa as the wavenumber or wavelength, respectively.

Introduction (Pharmaceutical substances)

Purpose and use of the spectra collection

Everyone knows how slim the chances are of recognizing and finding an unknown person purely on the basis of a few isolated details: Knowledge of his height, weight, age, the colour of his eyes and the style of his hair narrows the possibilities down slightly from a large number of individuals. It is the sum of numerous characteristics – most of which are registered unconsciously and can hardly be defined – that makes recognition possible. This is the reason why a clear picture – if possible together with exact details – is so much more informative.

The problems of spectroscopy are similar: Sometimes it is difficult to assess the UV or IR spectra of unknown substances from **theory**, and even the specialist can only rarely classify such substances to any degree of certainty.

Spectrophotometric assessment, both qualitative and quantitative, is based on comparison.

Absorption maxima and minima and specific UV absorption levels can be expressed in numerical data which provide valuable clues to recognition and will be used in the following tables, but it is also the width of an absorption band or a minimum range, or the steepness, the changes in direction or the minute shoulders of a curve that are often indispensable characteristics of a UV spectrum – just as the shape of the nose and mouth, curves and wrinkles really make up the face.

This being our experience, we have during the past few years compiled a collection of the UV spectra of important drugs, which means that we always have the most important comparative material at hand to make our search considerably easier. Subsequently the collection has been completed by the addition of the corresponding IR spectra.

Colleagues from other teams have been grateful for access to our collection. It thus seemed logical to reproduce a suitable selection of our collection and make it available to more users.

What is a »suitable selection«? There are lists of pharmaceutical agents of all kinds which name several thousand substances. There is a smaller number of substances with a broad connection with the therapeutic use of drugs; modern pharmacopoeias provide a guideline to these, and our choice is based on such reference books in general. Some of the substances listed in the pharmacopoeias have been omitted, while a greater number of new substances which do not appear in pharmacopoeias have been included because we felt that they were for various reasons important or interesting. It is probably inevitable that a selection of this nature will not meet with the unqualified approval of all its users. It is also quite obvious that this type of selection needs to be added to and updated constantly. We shall be pleased to receive suggestions for additions to the selection.

We have used non-registered recommended generic names (wherever possible, the INN) for the spectra, which we considered logical in view of the fact that the spectra are of active substances which have generic names, while preparations which take the form of pharmaceutical specialities have registered tradenames.

The ordinate for the UV spectra is the absorption value E_{1cm} of the solution measured (at a given concentration). We feel that this is the best presentation for the practical use of the spectra for the purpose of analysis, particularly for quantitative assessment. The commonly used ordinate log ϵ doubtless also has its advantages but we purposely chose to present the spectra \ast in their natural state \ast in the interest of practical applicability.

The IR spectra are presented as is common in percent transmission/wavelength or wave number using KBr disks 1.5 mg + 300 mg if not specified otherwise.

In our experience the collection can be of use in solving various problems:

It can contribute to the identification and classification of unknown substances. Large laboratories today obviously use NMR and mass spectroscopy to help them solve their identification problems, but it is frequently possible to ascertain the substance class of a sample with just a UV spectrum – if material and spectra are available for comparison. It is often possible to define the probable range of indications, as for example »analgesic«, »local anaesthetic«, »psychotropic« or similar. In that case, the field of major interest can be selected under »Indications«. The spectra sheets can be called up there. In some cases a pharmaceutical substance can be found in two or three different groups of indication.

Certainly for a reliable identification mostly UV absorption data are not sufficient. Therefore IR spectra must be included.

The UV absorption values in various solvents are expressed in E 1%/1cm and ϵ for quantification purposes. The values given in the table for each spectrum are those obtained in commercially available pure substance and checked on several photometers. Anyone who is familiar with this type of material will know that the values may differ slightly from those of other investigators.

Structural formulae are included to illustrate to the user the connections between the chemical structure and absorption properties so that he can develop his own spectrocognosia skills, which will allow him to recognize certain spectra or, in the case of unknown components, to classify a compound or recognize certain functional elements of the structure with some certainty, for example phenol hydroxy or aromatic amino groups, pyrazolones, xanthines, phenothiazines, etc. We have put together a few examples of this method of drug spectrocognition, showing well-defined structure/spectrum relationships. These examples are intended as a guide to practical application, but they are also intended as an incentive towards the user's development of his own **archives** for particular classes of substance.


As UV spectra are more distinct and readily classifiable we have used them firstly for classification.

Description and classification of substances is of course not always easy, and the information provided in a UV spectrum is not always sufficient for identification, but it should be possible with the aid of this type of picture and information to classify quite rapidly some substances not included in the selection from their UV spectra, according to the chemical (and sometimes also the therapeutic) group defined by their chromophores. For example, it is logical to assume an analgesic/antipyretic classification from a pyrazolone spectrum, while conversely if a substance is thought to be a local anaesthetic, the analyst will consider p-aminobenzoic acid esters or base-substituted anilides, and will read the spectrum from this aspect.

If there is no clue to the chemical group of a substance and the indication is not specified, the analyst can try to identify it by way of the list of absorption data, although it is not often possible to make a definitive identification of an unknown substance from the UV absorption data. It is, however, largely possible to classify material within this selection, together with reference to the spectra indicated by the specfinder list. For final identification the IR spectrum has to be included.

The first classification or recognition step (I) in our specfinder system is always the (clearly defined) UV absorption maximum with the longest wavelength in methanolic solution, which may also be the highest maximum of the whole range examined. It may in fact have a lower specific absorption level than a maximum with a shorter wavelength, but it is still considered as identification point I. If the list shows two or three maxima and one of those given under II is higher than I, the long-wave »principle maximum«, the higher one is underlined. It may sometimes be difficult to decide whether the maximum with the longest wavelength is »sufficiently clearly defined«, in which case the analyst is advised to look for a more clearly defined maximum at a shorter wavelength.

If a number of narrow bands in one group are relatively close, the most clearly defined band for this absorption range is the one that dictates identification. An example of this principle is tolbutamide, spectrum no. 508. Identification point I is the most clearly defined band in the long-wave range 263 nm (not 273 nm), and point II is the considerably more clearly defined short-wave band in methanol at $\underline{228}$ nm.

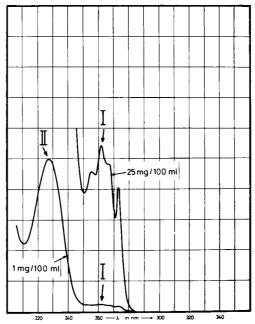

Tolbutamide / Attachment of bands

Another example of identification from the table: Clearly defined long-wave principle maximum I at 330 nm, second and third maxima at 254 and 219 = chloroquine phosphate (spectrum no. 908).

In some substances (e.g. particularly in substances with phenol OH or aromatic amino groups) aqueous alkaline or acid media cause the degree of dissociation to change, resulting in varying degrees of deformation of the spectrum compared with that in methanol. Changes of this nature are indicated by the following symbols in column IV of the specfinder table:

- ↑ increased absorption maximum↓ reduced absorption maximum
- $\uparrow\downarrow$ increase in one maximum, fall in another
- → bathochromic shift← hypsochromic shift
- increase + bathochromic shift
- ∠ fall + hypsochromic shift
- X considerable change in spectrum
- X 283 considerable change in spectrum + bearing a new maximum at 283 nm
- \rightarrow Z increasing degradation

Some substances only give clear spectra in diluted alkaline solution (0.1 M sodium hydroxide), and in such cases only that spectrum is shown. For those spectra a table has been included. Other substances are not sufficiently soluble in water, in which case only a methanol spectrum is shown.

Since there is naturally a certain degree of uncertainty in the expression of UV absorption maxima exactly to 1 nm, the user should in all cases take into consideration values up to \pm 2 nm of the value given.

Mostly quantitative absorption data (E 1%/1cm) are also required for classification or recognition; so of course solutions of known concentrations must be used as a basis.

With the exception of a few very similar compounds, this key permits perfectly serviceable allocation of the substances listed here. For reliable identification the IR spectra are available.

Introduction (Pharmaceutical substances) continued

Another important use of the collection is to supply absorption data on known drugs as basic information for analytical tasks, for example, to answer the following questions:

Is it possible to determine substance X by UV absorption methods, or is it promising under given conditions?

or

Is it possible to analyse substances X and Y (at given concentration ratios) by spectrophotometric methods either directly or with simultaneous measurement at 2 or 3 wavelengths or various pH values, or is preliminary separation (e. g. by chromatography) necessary or desirable?

These suggestions should suffice for the experienced, or even the new, analyst to find further uses for this collection of spectra.

We hope that the users of this book will find it helpful. Suggestions for changes or desirable additions coming from the analytical practice will be welcomed.

The technical assistance of Messrs. R. Merkel, M. Petry, and W. Weiss is gratefully acknowledged.

Further information

License conditions

Info

Publisher and publisher's products

Software license

ECV · Editio Cantor Verlag (ECV) software license

Please read the following license conditions carefully and thoroughly before using the software on your computer. By using the software or by opening the software packaging you explicitly agree to the following license conditions. If you do not agree to these license conditions, you must not use the software. In this case, you may immediately return the program package to ECV after purchasing or receiving it, and the purchase price will be refunded. The software is not sold, but licensed for the purpose of utilisation. It is merely the ownership of the storage medium (disk or CD), the manual and other related documents that passes to you.

1. Granting of a license

This license enables you

- to use a copy of the software on a single computer under the condition that the software is only used on a single computer at any one time. Using the software means that the software is either loaded into a computer's temporary memory (e. g. RAM) or a permanent memory (e. g. hard disk, CD-ROM). If you have acquired multiple software licenses, you may only use as many copies as licenses were acquired by you. You do not need additional licenses for a copy of the software installed on a generally accessible storage medium (e. g. a server). If the potential number of software users exceeds the number of acquired licenses, you must provide appropriate mechanisms or procedures to ensure that the number of persons using the software simultaneously does not exceed the number of licenses.
- to make a back-up copy of the software

2. Copyright

The software is protected by copyright. ECV claims the rights pursuant to the copyright. You are only authorized to use the software within the scope of permitted utilisation. Any other kind of utilisation is prohibited. The software contains copyrighted material as well as company secrets which you are obliged to keep. It is forbidden to decompile (§ 69 UrhG) or reassemble the software or convert it into any generally legible form and to modify, adapt, translate, hire out, lease out, lend, produce or publically display the software or parts thereof as well as products resulting herefrom.

The copyright especially refers to the program code, documentation, outer appearance, structure and organisation of program files, the program name, logos and other forms of presentation within the software. Without the explicit and written consent of ECV, copying, processing, any other kind of conversion, modification or reproduction of the contents of the software is prohibited.

The manual as well as other documents related to the software are protected by copyright. It is forbidden to copy, process, convert, modify or distribute the printed material. Any offenses will be prosecuted according to civil or criminal law.

You must not hire out, lease out, sell sub-licenses, lend the software or parts of it or agree to copying the software onto another user's computer. However, you may transfer all your rights concerning the usage of the software to another natural or legal person under the condition that a) you transfer to that person this agreement, the software and any other software or hardware delivered with the software or preinstalled on it, including any copies, updates and previous versions as well as any copies of the documents converted into any other formats as well as the entire documentation (manuals etc.), b) you do not retain any copies, including back-up copies or any other copies stored on a computer and c) the recipient approves of the conditions stipulated in this agreement or any other conditions under which you have legally obtained the software license. Training, preliminary or sample copies of the software must never be transferred.

3. License duration

The license is granted for an unlimited period of time. The license automatically loses effect without prior notice if you offend against any term of this agreement. On terminating the agreement, you are obliged to destroy the software and any copies thereof. You may terminate the license agreement at any time by destroying the software including all copies.

4. Deficiencies, warranty

4.1. The parties are aware that software cannot be produced completely faultlessly.

Information contained in the manual/documentation and/or in the advertising material relating to a product's upgradability or available accessories is not binding, as the products are subject to continuous modification, and information may also pertain to future developments.

Manuals and documentation other than the documents/program description supplied with the software and user guidance implemented in the software and/or online help shall only be provided if both parties have expressly agreed in writing. In case of such an explicit agreement, no requirements are made concerning the contents, language and volume of a manual and/or documentation to be supplied explicitly, and supply of a brief description shall be considered sufficient, unless otherwise agreed in writing by both parties.

- 4.2. No liability is assumed by ECV for the licensed software achieving certain performance results or being able to be utilised for the intended purpose or for compatibility with the software already employed by the user. Especially the risk of economic usability exclusively lies with you.
- 4.3. Any products or services delivered must be checked immediately upon receipt. Any claims must be made in writing without any delay. Complaints about obvious faults may only be filed within a period of two weeks.

When making a claim in writing, precise information must be given on the content and to which end the software should be utilized according to the contract, which and how many steps have been taken and, if available, which error messages were output by the software.

- 4.4. With the exception of the purchasing of consumer goods pertaining to new goods, to which the legal two-year period of limitation is applicable, the right to make a claim with regard to a deficiency expires after one year from the beginning of the legal period of limitation. The legal periods of limitation apply in case of a violation of life, body or health or any other damage resulting from intentional or grossly negligent breach of duty or in case of malicious silence with regard to a deficiency for which ECV must be made responsible. In case of a deficiency consisting in the right in rem of a third party to be surrendered the purchased item (§ 438 I 1 BGB), the period of limitation is five years.
- 4.5. In other cases, the legal regulations regarding the warranty for deficiencies apply.
- 4.6. However, when entrepreneurs are concerned, ECV reserves the right to choose whether to eliminate a fault or deliver a faultless item (§ 439 I BGB) in case of justified claims for the rectification of faults.
- 4.7. If ECV is liable, the following applies: Liability for less serious negligence is only assumed in case of delay, impossibility and neglect of an essential duty stipulated in the agreement. In such cases, liability is furthermore restricted to the predictable typical damage (typical average damage as stipulated in the agreement).

If the licensee is an entrepreneur, ECV assumes no liability either for indirect damage, consequential damage, lost profits or any other kind of economic loss.

This limited warranty does not hold true for damage resulting from the violation of life, body or health. In other cases, the legal liability of ECV remains unaffected, especially the liability for malicious aforethought, gross negligance, malicious silence with regard to a deficiency, any guarantee assumed and in accordance with the product liability law.

Software-License continued

- 4.8. This limited warranty also holds true for third parties considered in the scope of protection of the contractual relation as well as in favour of legal representatives, employees, collaborators or any other ECV helpers or assistants.
- 4.9. If ECV is required to assume liability or warranty claims are filed, it must be considered that the user may also be responsible to a certain extent, especially in case of inadequate error reports or data storage.

5. Other

This license agreement is subject to the law of the Federal Republic of Germany. If terms of this license agreement are or become ineffectual in full or in part, this shall have no effect on other terms. A term rendered ineffectual must rather be replaced by such a term the object and purpose of which is as similar as possible to that of the ineffectual term. No other agreements have been made. Modifications to this license agreement must be submitted in writing. This also applies to cancellation of the written-form requirement.

ECV · Editio Cantor Verlag für Medizin und Naturwissenschaften GmbH, Aulendorf (Germany)

UV and IR Spectra of Pharmaceutical Substances and IR Spectra of Pharmaceutical and Cosmetic Excipients

H.-W. Dibbern, R. M. Müller, E. Wirbitzki (Eds.)



© 2002

Editio Cantor Verlag für Medizin und Naturwissenschaften GmbH

P.O. Box 1255 · 88322 Aulendorf (Germany)
Phone +49 (0)7525-940 0 · Fax +49 (0)7525-940 180 · e-mail info@ecv.de · http://www.ecv.de/

Important notes

- The absence of the registered trademark symbol ® after a name does not mean that the particular name is not protected by trademark.
- The editors of this work have made every effort to provide correct data. The measuring instruments were state of the art at the time they were used to compile the respective data. They were calibrated at regular time intervals according to the manufacturers' instructions and in compliance with pharmacopoeial specifications. The pharmaceutical substances used for compilation of the data complied with the pharmacopoeial specifications resp. the quality requirements set out by the manufacturers. Regardless of taking greatest care neither the editors nor the publisher can assume any guaranty for the data.
- The work and any part of it are protected by copyright. Any exploitation, especially any reproduction other than that permitted by law, is strictly prohibited.
- Adobe, Acrobat, PageMaker, Type Manager, Illustrator, and PostScript are trademarks of Adobe Systems Incorporated. Parts of Adobe File Utilities were licensed by Adobe Systems Incorporated. Copyright © 1986 – 1997. All rights reserved.
- All other trade names and product names, typefaces, company names and logos of companies are trademarks or registered trademarks of the respective companies.

Design, processing, conception, and realization by TD-PIRSCH 2002

Info: http://www.tdpirsch.de/



Further information [Publisher and publisher's products]

Please refer to our homepage http://www.ecv.de/ for information about the publisher and the publisher's products.

Weitere Informationen

Lizenzbedingungen

Info

Verlag und Programm

Software-Lizenz

Software-Lizenz des ECV · Editio Cantor Verlages (ECV)

Lesen Sie nachfolgende Lizenzbedingungen aufmerksam und sorgfältig durch, bevor Sie die Software auf Ihrem Computer einsetzen. Durch Verwendung der Software bzw. durch Öffnen der Software-Verpackung erklären Sie Ihr ausdrückliches Einverständnis mit den nachstehenden Lizenzbestimmungen. Für den Fall, daß Sie mit diesen Lizenzbedingungen nicht einverstanden sind, dürfen Sie die Software nicht verwenden. In diesem Fall können Sie das Programmpaket unverzüglich nach Erwerb oder Erhalt an ECV zurücksenden und erhalten den Kaufpreis rückerstattet. Die Software wird nicht verkauft, sondern lizenziert zum Zwecke der Nutzung. Eigentum erhalten Sie nur am Speichermedium (Diskette oder CD) sowie am Handbuch sowie den sonstigen zugehörigen Schriftdokumenten.

1. Einräumung einer Lizenz

Diese Lizenz erlaubt Ihnen:

- die Benutzung einer Kopie der Software auf einem Einzelcomputer unter der Voraussetzung, daß die Software zu jeder Zeit auf nur einem einzigen Computer verwendet wird. Die Benutzung der Software bedeutet, daß die Software entweder in einem temporären Speicher (z. B. RAM) eines Computers oder auf einem permanenten Speicher (z. B. Festplatte, CD-ROM) geladen ist. Wenn Sie Mehrfachlizenzen für die Software erworben haben, dürfen Sie immer nur höchstens so viele Kopien in Benutzung haben, wie Lizenzen von Ihnen erworben wurden. Sie benötigen keine zusätzliche Lizenz für eine Kopie der Software, die auf einem allgemein zugänglichen Speichermedium (z. B. Server) selbst installiert ist. Wenn die voraussichtliche Zahl der Benutzer der Software die Zahl der erworbenen Lizenzen übersteigt, so müssen Sie angemessene Mechanismen oder Verfahren bereithalten, um sicherzustellen, daß die Zahl der Personen, die die Software gleichzeitig benutzen, nicht die Zahl der Lizenzen übersteigt.
- eine Kopie der Software für Sicherungszwecke herzustellen.

2. Urheberrecht

Die Software ist urheberrechtlich geschützt. Die aus dem Urheberrecht resultierenden Rechte stehen ECV zu. Sie sind nur berechtigt, die Software im Rahmen der erlaubten Nutzung zu verwenden. Jede darüber hinausgehende Nutzung ist untersagt. Die Software enthält urheberrechtlich geschütztes Material sowie Betriebsgeheimnisse, zu deren Wahrung Sie sich verpflichten. Es ist verboten, die Software zu dekompilieren (§ 69 UrhG), rückassemblieren oder auf andere Weise in allgemein lesbare Form umzuwandeln sowie Software oder Teile der Software sowie hieraus abgeleitete Produkte zu ändern, anzupassen, zu übersetzen, zu vermieten, zu verleasen, zu verleihen, herzustellen oder öffentlich wiederzugeben.

Das Urheberrecht umfaßt insbesondere den Programmcode, die Dokumentation, das Erscheinungsbild, die Struktur und Organisation der Programmdateien, den Programmnamen, Logos und andere Darstellungsformen innerhalb der Software. Jede Vervielfältigung, Nutzung, Bearbeitung, sonstige Umgestaltung, Änderung oder Wiedergabe des Inhaltes der Software ohne ausdrückliche schriftliche Einwilligung von ECV ist untersagt.

Das Handbuch sowie sonstige zur Software gehörende Schriftstücke sind urheberrechtlich geschützt. Jede Vervielfältigung, Bearbeitung, sonstige Umgestaltung, Änderung oder öffentliche Wiedergabe des Schriftmaterials ist verboten und wird zivil- und strafrechtlich verfolgt.

Sie dürfen die Software weder in Teilen noch als Ganzes vermieten, verpachten, unterlizenzieren, verleihen oder dem Kopieren der Software auf den Computer eines anderen Benutzers zustimmen. Sie dürfen jedoch alle Ihre Rechte zur Verwendung der Software an eine andere natürliche oder juristische Person unter der Voraussetzung übertragen, daß a) Sie den vorliegenden Vertrag, die Software und sonstige Software oder Hardware, die mit der Software geliefert oder auf dieser vorinstalliert ist, einschließlich aller Kopien, Updates und früherer Versionen sowie aller Kopien der Schriftstücke, die in andere Formate konvertiert wurden, sowie die gesamte Dokumentation (Handbücher etc.) an diese Person übertragen, b) Sie keine Kopien, einschließlich Sicherheitskopien oder sonstigen Kopien, die auf einem Computer gespeichert sind, zurückbehalten und c) der Empfänger in die Bestimmungen dieses Vertrages sowie sonstige Bestimmungen eingewilligt hat, nach denen Sie die Softwarelizenz rechtmäßig erworben haben. Schulungs-, Vorab- oder Musterkopien der Software dürfen niemals übertragen werden.

3. Dauer der Lizenz

Die Einräumung der Lizenz erfolgt zeitlich unbefristet. Die Lizenz verliert automatisch ihre Wirksamkeit, ohne daß es einer Kündigung bedarf, wenn Sie gegen irgendeine Bestimmung dieses Vertrages verstoßen. Im Falle der Beendigung sind Sie verpflichtet, die Software sowie alle Kopien der Software zu vernichten. Sie können den Lizenzvertrag jederzeit dadurch beenden, daß Sie die Software einschließlich aller Kopien vernichten.

4. Mängelrechte, Haftung

4.1. Die Parteien sind sich darüber einig, daß Software nicht völlig fehlerfrei erstellt werden kann.

Angaben im Handbuch/in der Dokumentation und/oder im Werbematerial, die sich auf Erweiterungsmöglichkeiten eines Produkts beziehen oder auf verfügbares Zubehör, sind unverbindlich, insbesondere weil die Produkte ständiger Anpassung unterliegen und sich die Angaben auch auf zukünftige Entwicklungen beziehen können.

Die Lieferung von Handbüchern und Dokumentationen über das mit der Software ausgelieferte Schriftmaterial/Programmbeschreibung und die in die Software implementierte Benutzerführung und/oder Online-Hilfe hinaus wird nur dann geschuldet, wenn dies ausdrücklich schriftlich zwischen den Parteien vereinbart worden ist. Im Falle einer solchen ausdrücklichen Vereinbarung sind Anforderungen hinsichtlich Inhalt, Sprache und Umfang eines ausdrücklich zu liefernden Handbuches und/oder Dokumentation nicht getroffen, und die Lieferung einer Kurzanleitung ist ausreichend – es sei denn, daß die Parteien schriftlich weitere Spezifikationen vereinbart haben.

- 4.2. ECV haftet nicht dafür, daß die lizenzierte Software bestimmte Leistungsergebnisse erzielt oder für den von Ihnen beabsichtigten Zweck eingesetzt werden kann oder mit beim Anwender vorhandener Software zusammenarbeitet. Insbesondere liegt das Risiko der wirtschaftlichen Verwertbarkeit allein bei Ihnen.
- 4.3. Alle Lieferungen und Leistungen sind sofort nach Empfang zu prüfen. Beanstandungen sind unverzüglich in Textform mitzuteilen. Offensichtliche Mängel können nur binnen einer Frist von zwei Wochen geltend gemacht werden.

Im Rahmen der schriftlichen Mängelrüge sind konkrete Angaben dahingehend zu machen, mit welchem Inhalt und Ziel die Software vertragsgemäß betrieben werden sollte, welche und wie viele Arbeitsschritte vorgenommen worden sind und, soweit vorhanden, mit welchen Fehlermeldungen die Software reagiert hat.

- 4.4. Mit Ausnahme beim Verbrauchsgüterkauf über neue Sachen, für den die gesetzliche zweijährige Verjährungsfrist gilt, verjähren Mängelrechte nach Ablauf eines Jahres ab dem gesetzlichen Verjährungsbeginn. Bei einer von ECV zu vertretenden Verletzung von Leben, Körper oder Gesundheit sowie bei einem sonstigen Schaden, der auf einer vorsätzlichen oder grob fahrlässigen Pflichtverletzung beruht sowie bei arglistigem Verschweigen eines Mangels gelten die gesetzlichen Verjährungsfristen. Liegt ein Mangel vor, der in einem dinglichen Recht eines Dritten auf Herausgabe der Kaufsache besteht (§ 438 I 1 BGB), beträgt die Verjährungsfrist fünf Jahre.
- 4.5. Im übrigen gelten für die Mängelrechte die gesetzlichen Regelungen.
- 4.6. Gegenüber Unternehmern steht bei berechtigten Nacherfüllungsansprüchen das Wahlrecht, den Mangel zu beseitigen oder eine mangelfreie Sache zu liefern (§ 439 I BGB), jedoch ECV zu.
- 4.7. Haftet ECV, so gilt folgendes: Die Haftung bei einfacher Fahrlässigkeit besteht nur bei Verzug, Unmöglichkeit und Verletzung einer wesentlichen Vertragspflicht. In diesen Fällen ist die Haftung ferner auf den voraussehbaren typischen Schaden (vertragstypischer Durchschnittsschaden) begrenzt.

Ist der Lizenznehmer ein Unternehmer, haftet ECV ferner nicht für mittelbare Schäden, Mangelfolgeschäden, entgangenen Gewinn und sonstige Vermögensschäden.

Diese Haftungsbeschränkung gilt nicht für Schäden aus der Verletzung des Lebens, des Körpers oder der Gesundheit. Im übrigen bleibt die gesetzliche Haftung des ECV unberührt, insbesondere die Haftung für Vorsatz, grobe Fahrlässigkeit, arglistiges Verschweigen eines Mangels, eine etwa übernommene Garantie sowie nach dem Produkthaftungsgesetz.

- 4.8. Die Haftungsbeschränkung gilt auch gegenüber Dritten, die in den Schutzbereich der Vertragsbeziehung einbezogen werden sowie zugunsten der gesetzlichen Vertreter, Arbeitnehmer, Mitarbeiter und sonstigen Verrichtungs- und Erfüllungsgehilfen des ECV.
- 4.9. Im Falle einer Inanspruchnahme von ECV aus Gewährleistung oder Haftung ist ein Mitverschulden des Anwenders angemessen zu berücksichtigen, insbesondere bei unzureichenden Fehlermeldungen oder unzureichender Datensicherung.

5. Sonstiges

Dieser Lizenzvertrag unterliegt dem Recht der Bundesrepublik Deutschland. Für den Fall, daß Bestimmungen dieses Lizenzvertrages ganz oder teilweise unwirksam sind oder werden, so berührt dies die Wirksamkeit der übrigen Bestimmungen nicht. Die unwirksame Bestimmung ist vielmehr durch eine solche zu ersetzen, die dem Sinn und Zweck der unwirksamen Bestimmung möglichst nahekommt. Nebenabreden sind nicht getroffen. Änderungen dieser Lizenzvereinbarung bedürfen der Schriftform. Gleiches gilt für die Aufhebung dieser Schriftformklausel.

ECV · Editio Cantor Verlag für Medizin und Naturwissenschaften GmbH, Aulendorf

UV and IR Spectra of Pharmaceutical Substances and IR Spectra of Pharmaceutical and Cosmetic Excipients

H.-W. Dibbern, R. M. Müller, E. Wirbitzki (Hrsg.)



© 2002

Editio Cantor Verlag für Medizin und Naturwissenschaften GmbH

Postfach 1255 · 88322 Aulendorf (Germany)

Telefon +49 (0)7525-940 0 · Fax +49 (0)7525-940 180 · e-mail info@ecv.de · http://www.ecv.de/

Wichtige Hinweise

- Das Fehlen des Symbols ® nach Namen bedeutet nicht, daß der Name nicht durch Warenzeichen geschützt ist.
- Die Herausgeber dieses Werkes haben sich intensiv bemüht, korrekte Angaben zu machen. Die für die Erstellung der Daten verwendeten Meßgeräte entsprachen zu jeder Zeit dem aktuellen Stand der Technik. Die Kalibrierung erfolgte regelmäßig nach den Anweisungen der Gerätehersteller sowie nach den geltenden Vorschriften der Arzneibücher. Die für die Erstellung der Daten eingesetzten Wirkstoffe entsprachen den Arzneibuchspezifikationen bzw. den Qualitätsanforderungen der Hersteller. Ungeachtet größtmöglicher Sorgfalt können Herausgeber und Verlag für die Daten keine Gewähr übernehmen.
- Das Werk und seine Teile sind urheberrechtlich geschützt. Jede Verwertung insbesondere Vervielfältigung in anderen als den gesetzlich zugelassenen Fällen ist strengstens untersagt.
- Adobe, Acrobat, PageMaker, Type Manager, Illustrator, PostScript sind Warenzeichen von Adobe Systems Incorporated. Teile der Adobe File Utilities wurden von Adobe Systems Incorporated lizensiert. Copyright © 1986–1997. Alle Rechte vorbehalten.
- Alle anderen Handels- und Produktbezeichnungen, Schriften sowie Firmennamen und -logos sind Warenzeichen oder eingetragene Warenzeichen der jeweiligen Unternehmen.

Entwurf, Aufbereitung, Konzeption und Realisation by TD-PIRSCH 2002

Info: http://www.tdpirsch.de/



Weitere Informationen [Verlag und Programm]

Informationen über den Verlag und sein Programm finden Sie auf unserer Homepage: http://www.ecv.de/

Spectra

Pharmaceutical Substances

Pharmaceutical and Cosmetic Excipients

Pharmaceutical Substances

- On the left part of the screen, you can select the active pharmaceutical ingredient or excipient category (in the case of the active pharmaceutical ingredients, a letter of the alphabet is also displayed) and open the desired spectra sheet by clicking on it.
- Sie können in der Leiste links die Rubrik Wirkstoff oder Hilfsstoff auswählen (bei den Wirkstoffen zusätzlich einen Buchstaben aus dem Alphabet) und anschließend das gewünschte Spektrenblatt durch Anklicken öffnen.

Pharmaceutical and Cosmetic Excipients

- On the left part of the screen, you can select the active pharmaceutical ingredient or excipient category (in the case of the active pharmaceutical ingredients, a letter of the alphabet is also displayed) and open the desired spectra sheet by clicking on it.
- Sie können in der Leiste links die Rubrik Wirkstoff oder Hilfsstoff auswählen (bei den Wirkstoffen zusätzlich einen Buchstaben aus dem Alphabet) und anschließend das gewünschte Spektrenblatt durch Anklicken öffnen.

Indications

- On the left part of the screen, you can select a therapeutic indication and open one of the displayed spectra sheets by clicking on it.
- Sie können in der Leiste links eine Indikation auswählen und anschließend eines der aufgeführten Spektrenblätter durch Anklicken öffnen.

UV Spectrocognosia

The following pages do not attempt a general explanation of spectra or a study of all possible structural effects; nor do they set out to describe the characteristics of all classes of substances. The explanations on UV spectra hereafter are merely intended to teach non-specialist users of spectroscopy techniques to recognize certain simple types of spectra which occur frequently in drug analysis so that they will be capable of at least limited identification of the class of substance in simple cases. When evaluating unknown spectra it must always be borne in mind that complex structural factors sometimes make otherwise clear, typical bands difficult to recognize or even unrecognizable. Knowledge of the major basic types is nevertheless valuable for those who analyse drugs.

It is intended to give the user an introduction to UV spectrocognosia of pharmaceuticals and to teach him how to apply his knowledge of such »family characteristics«. Some users will mainly be concerned with specific classes of substance, and will develop their own spectrocognosia. The following examples may be of assistance.

Phenyl groups without other chromophoric groups

This type of spectrum occurs frequently in the pharmaceutical field. The chromophore has 3–5 relatively close bands at 246–268 nm. Clear distinction between the narrow bands depends not only on the quality of the photometer but also on the substance in hand. Good differentiation depends on the analysis of the sample at the optimum concentration. The wavelengths of the bands can be seen in figure 1.

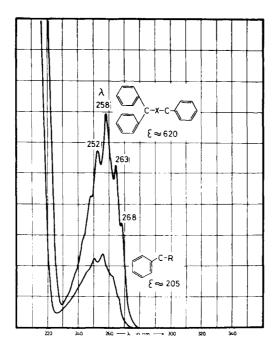


Fig. 1: Phenyl- without functional groups.

The specific absorption of the pure »phenyl chromophore« is quite low: for \bigcirc - R type compounds $\epsilon \approx 200$ (in methanol). If the molecule has 2 or 3 phenyl radicals, ϵ increases to about 400 or 600. Thus the E 1%/1cm for many drugs is around 5–20, i. e. measurements should be taken in concentrations of about 20–50 mg/100 ml.

The Specfinder Table contains a large number of examples at about 257 nm. Phenmetrazine hydrochloride [007] has a slight hypsochromic shift up to 255 nm.

In the case of benzyl penicillin [402] the short wavelength part of the spectrum is »blurred«, but the narrow bands at 263 and 257 are still well defined.

Phenol OH groups

	Spectrum	Absorpt.	Shoulder	Change in
	no.	max.	at λ max.	alkaline
		in methanol	+ 8nm	solution
Phenol	1210	273 nm	+	7
Tyramine	2709	277 nm	+	7
Salicylic acid	1105	302 nm	Ø	Ľ
Salicylic acid methylester	703	305 nm	Ø	7
Salicylamide	118	288 nm	Ø	7
p-Hydroxybenzoic acid	1205	252 nm	Ø	7
p-Hydroxybenzoic acid ester	1207	256 nm	Ø	7
Paracetamol	112	247 nm	Ø	И

Phenol OH groups have a very characteristic appearance in the UV spectrum. If there are no other functional groups on the benzene ring, the absorption maximum in methanol is at about 274 ± 3 nm, with a very typical »shoulder« at 282 ± 3 nm (see figure 2). This applies not only to phenol [1210], thymol [1212] and resorcinol [1104] but also to many antihypotensive agents and sympathicomimetic drugs from Group 27. The phenol oestrogens have a bathochromic shift of about 8 nm (estradiol [2202], estradiol 17β -cipionate [2203], estrone [2204], ethinylestradiol [2205]). Ortho-diphenols (corbadrin [2701], etc.) do not have the shoulder, and neither do all compounds which have other functional groups on the benzene ring. There are always considerable shifts towards long or short wavelengths with bifunctional compounds of this type.

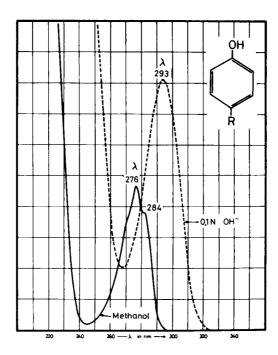


Fig. 2: Phenol. OH- without other functional groups.

Almost all phenol compounds in aqueous alkaline solution display a marked bathochromic shift and usually (but not always) an increase in absorption (symbol \nearrow) comparison with solution in methanol, as result of ionisation. Salicylic acid [1105] is an exception: the influence of the free ortho-carboxyl group predominates. Salicylic acid esters [703] display the alkali shift again. Ortho-diphenols are unstable in alkaline solution and can therefore hardly be measured.

Note: Bathochromic shifts \rightarrow and \nearrow in alkaline solution are not exclusive to phenols, but also occur in hydroxypyridines [2806] and ketones capable of forming enols: pyridone [909] and benzoxazolinone [2402], in benzodiazepinones [2604] and also in nitro-compounds (metronidazole [916], nifurprazine [917], nitrofural [918], nitrofurantoin [919]).

UV Spectrocognosia continued

Nevertheless, the shoulder at λ max + 8 nm plus the 7 shift can be regarded as good indication of phenol OH, particularly when there are no other functional groups on the ring. A 7 shift on its own provides only a limited indication of identity.

The intensity of the phenol band at about 275 nm is moderately pronounced: ϵ is of the order of 2000 (depending on other substituents), and in the case of ortho-diphenols it is about 3000.

Phenol ethers

Phenol ethers in methanol, particularly those without any other functional groups, have a relatively similar UV spectrum to phenols: λ max at about 270–285 nm, sometimes with no marked shoulder (meticillin sodium salt [408], verapamil hydrochloride [1416]), sometimes with a faintly defined shoulder (guaifenesin [2406]) and sometimes with a more marked narrow band at λ max + approx. 8 nm (chlorphenesin carbamate [1101]).

There is of course no bathochromic shift here in alkaline solutions (difference between morphine hydrochloride [111] and codeine phosphate [801]). In many compounds the relatively weak phenol ether band is considerably overshadowed and barely discernible (mepyramine maleate [306], noscapine hydrochloride [805], papaverine hydrochloride [2408]).

The intensity of phenol ether absorption is lower than that of the phenols: ε is around 1000.

Aromatic carboxylic acids

Aromatic carboxylic acids (benzoic acid [2902], aspirin [102]) have a moderately characteristic chromophore: the absorption maxima are around 273 ± 3 nm and 227 ± 3 nm and can only be seen as an indication of identity together. Ionisation in NaOH reduces the absorption compared with the solution in methanol. Other chromophores present at the same point usually cover the aromatic carboxyl absorption.

Among the pyridine carboxylic acids, the carboxyl function is only more marked in the spectrum of isonicotinic acid [2907].

Aromatic amines

Aromatic amines are among the strongest and most important chromophores in drug analysis. They have an absorption maximum of around 290 nm, but this shifts considerably in some cases according to substituent and functional group (see figure 3). The high intensity (in methanol) has a characteristic $\epsilon \approx 20000$, and this chromophore therefore shows through very clearly. Another extremely marked characteristic is the suppression of this band almost to zero, as the H⁺-concentration in the solution increases, i. e., as ionisation increases (it is possible to determine the pKa of aromatic amines on the basis of this). As long as other molecular components do not interfere at short wavelengths, a band appears or becomes more prominent at about 220–240 nm as ionisation increases (benzocaine [1701], tetracaine hydrochloride [1709], acriflavinium chloride [902], sulfaguanidine [928]). In 0.1 M sodium hydroxide solution the main band generally shifts slightly towards lower wavelengths.

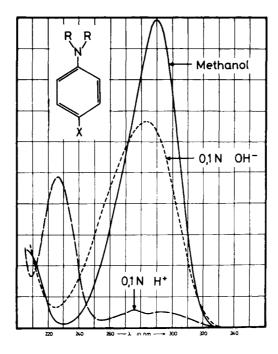


Fig. 3: Aromatic amines.

In principle, primary aromatic amines have the same type of absorption characteristics as secondary and tertiary aromatic amines (with aliphatic radicals).

Sulfonamides

Sulfonamides (phthalylsulfathiazole [920], sulfacetamide [921], sulfaquinoxaline [922], etc.) with chemotherapeutic activity are without exception 4-aminobenzene sulfonamides (sulfanilamides). Their UV absorption is thus mainly determined by the aromatic amino group. Sulfanilamides have one major band, usually at about 270 nm.

This is a considerably shorter wavelength than seen in 4-aminobenzoic acid esters, for example. The fall in absorption in acid solution is not equally pronounced in all cases.

Naturally the position of the major band and the ionisation change are affected by the substituents in the SO_2 -NH group and amino group: Strong chromophoric components like phthalyl- and quinoxalyl- [920, 922] show through clearly as a result of this. This means that the UV spectra of chemotherapeutic sulfonamides are not always immediately recognizable as typical.

Sulfonamide diuretics

Sulfonamide diuretics (cyclothiazide [1307], furosemide [1309], hydrochlorothiazide [1310], bendroflumethiazide [1315], hydroflumethiazide [1316], butizide [1317], cyclopenthiazide [1318]) with the general formula

are readily recognizable from a strong absorption band at 270–278 nm ($\epsilon \approx 20000$), a weaker band at 315–325 nm ($\epsilon \approx 3000$ –4000) and a low minimum at 240 nm (see figure 4). The longwave band in furosemide [1309] undergoes a further bathochromic shift (\rightarrow 336 nm) under the influence of the furyl substituent.

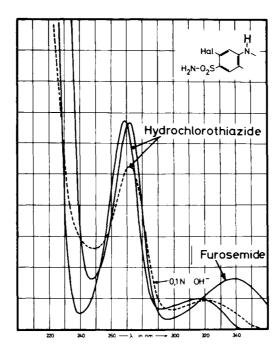


Fig. 4: Sulfonamide diuretics.

Chlorothiazide [1305] (without N-**H** on the benzene sulfonamide ring), bumetanide [1302] and chlortalidone [1303] do not belong to the above type of spectra.

Anilides

For the purpose of drug spectrocognition, acid amides of aromatic amines have to be subdivided into two largely different groups:

- A) Simple acid anilides of the acetanilide/phenacetine type [101, 114] and
- B) anilides substituted at 2.6 in the aromatic ring, of the lidocaine/butanilicaine type [1707, 1703].

Type A substances have a strong absorption band at 240–250 nm ($\epsilon \approx 1500$), far stronger in methanol than in aqueous (alkaline) solution or acid (see figure 5). Paracetamol [112] belongs to this group; it can be distinguished from phenacetin [114] easily because of a $\mbox{\em s}$ shift in 0.1 M OH $^-$ (due to ionization of the phenol group).

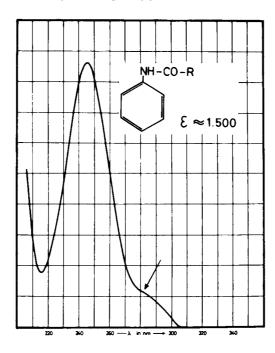


Fig. 5: Type A acetanilides.

Phenacetin [114], bucetin [104] and paracetamol [112] have a relatively characteristic shoulder at 280–290 nm.

NB: Acetohexamide [501] has a very similar spectrum despite the fact that its structure is entirely different.

2,6 substituted anilides of type B have only a weak band at about 265 nm ($\epsilon \approx 350$) with a narrow band at λ max + 8 nm (see figure 6).

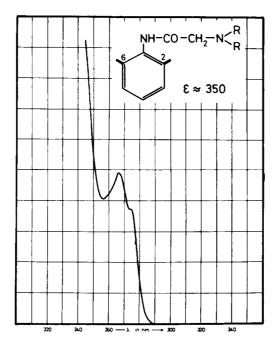


Fig. 6: Type B acetanilides.

The absorption in aqueous solution is slightly higher (\uparrow) than in methanol. Around 245 nm, type B compounds only show non-specifically increasing absorption with decreasing UV wavelength but no marked band as described for type A.

Pyrazolone compounds

Pyrazolone compounds with pharmaceutical applicability are, in UV spectrophotometry, a relatively limited, clearly defined group in which phenylbutazone [701] can be included. They are characterised by 2 relatively broad maxima with different wavelengths and intensities, 234–246 nm and 265–275 nm (see figure 7). In the base-substituted compounds (dimethylamino phenazone [103], metamizole sodium [109]) in acid solution the longwave maximum occurs at approximately the point of the minimum in methanol solution. The curves are slightly flattened in alkaline solution. ϵ is around 10000 in methanolic solution.

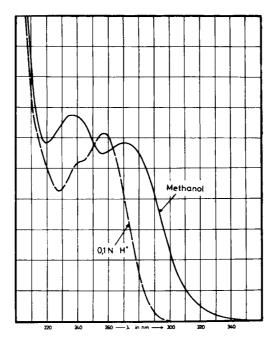


Fig. 7: Pyrazolones.

Pyridyl

The absorption of the pyridyl radical is considerably higher than that of unsubstituted phenyl. There is a clear band at about 261 nm with smaller bands (shoulders) at 269 and 265 nm (see figure 8).

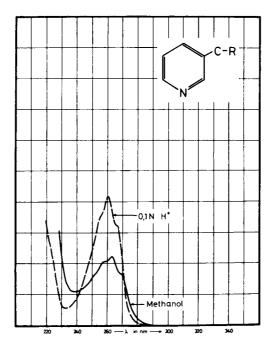


Fig. 8: Pyridyl-.

An important characteristic is the considerable increase in the band in 0.1 M acid (for example pheniramine maleate [308], nicotine [3003]). Some of these characteristics persist in the presence of functional groups, as in nicotinic acid and nicotinic acid amide [1802, 2805]. The narrow bands disappear in isonicotinic acid [2907]. In hydroxy-pyridine compounds [2806] the character of the spectrum changes completely, similar in some respects to the phenols (e. g. \rightarrow shift in alkaline solution).

Indoles

Indoles have a very characteristic spectrum which is readily recognized: it has a peak at 280 nm, a shoulder at about 274 nm and a clear narrow band at 289 ± 1 nm (see figure 9). Perfect examples of this can be seen in tryptamine [2909] and tryptophan [3006].

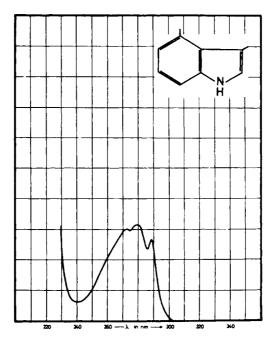


Fig. 9: Indolyl.

These indole characteristics also occur in hydrated ergot-alkaloids but not in lysergic acid derivatives.

The spectrum for pindolol [2008] still shows the narrow band at 288 nm very clearly, while the indole peak at 280 nm is only a shoulder (the pindolol maximum at 265 nm is also affected by the indole ether group).

These examples show that the analyst must not only look at the maximum when attempting to classify a compound.

Xanthines

The xanthines used in medicine are quite a homogenous group as regards UV spectra. All the compounds listed (caffeine [002], pentoxifylline [006], diprophylline [1406], etofylline [1407], theobromine [1414], theophylline [1415]) have a good, constant absorption band at 272 \pm 2 nm; ϵ approximately 9000–10000 in methanol; minimum 242 \pm 2 nm (see figure 10).

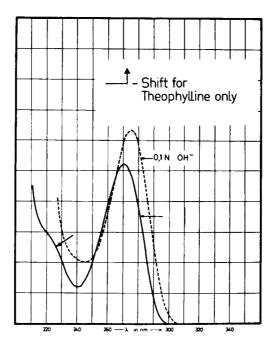


Fig. 10: Xanthines.

Xanthines can be identified with certainty because they have a weak but very typical shoulder at 225-230 nm. Only theophylline [1415], which is the only notable ionisating compound in this series, exhibits a relatively pronounced \nearrow shift in alkaline solution.

Phenothiazines

Phenothiazines can be readily recognized from the UV spectrum: they have a characteristic, broader, lower band at long wavelengths, a clear band at about 250 nm (ϵ approx. 30000), and characteristic minima (see figure 11).

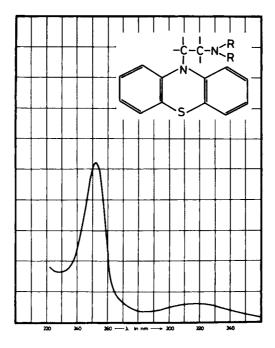


Fig. 11: Phenothiazines.

Neuroleptics with basic side-chains (chlorpromazine hydrochloride [2102], profenamine hydrochloride [2107]) shift slightly towards shorter wavebands than the parent substance (phenothiazine [203]).

Benzodiazepines

Benzodiazepines of the familiar 1,4-benzodiazepin-2-one series have in acid solution a characteristic, rather asymmetrical absorption band at about 283 nm ($\epsilon \approx 12000$) and a flatter maximum at 360 nm (see figure 12). Oxazepam [2605] can be distinguished from diazepam [2604] by the fact that in alkaline solution there is a definite shift of the maximum to a higher wavelength (in methanol 315, in alkali \rightarrow 344 nm).

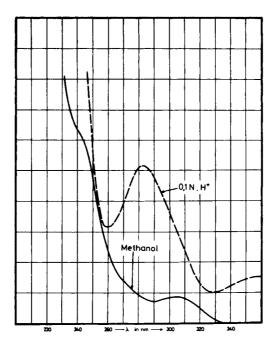


Fig. 12: 1,4-Benzodiazepine-2-ones.

The spectra of chlordiazepoxide hydrochloride [2602] and clobazam [2603] (1,5-benzodiazepine-2,4-dione) are markedly different from those of the series just described.

Barbiturates

UV spectrophotometry of the barbiturates is not very informative. UV spectra with any degree of clarity can only be obtained in alkaline solution at rather short wavelengths. The extended chromophore in phenobarbital [1508] has a maximum at 252 nm, as compared with the λ max. of about 242 nm for purely aliphatic substituted substances (see figure 13).

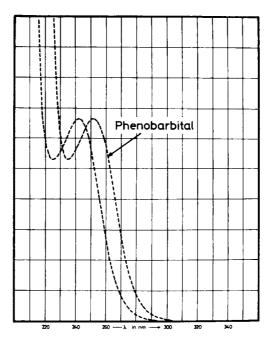


Fig. 13: Barbiturates in 0.1 M OH⁻.

Steroid hormones

Steroid hormones fall readily into 3 categories according to their UV spectrum:

1. The estrogens (estradiol 17β -cipionate [2203], estrone [2204], ethinylestradiol [2205]) have the absorption pattern already described for the phenols because they have a phenol OH group. The only difference is that the spectra of the estrogens undergo a slight shift towards longer wavelengths (λ max. approx. 280 nm) than the simple phenols (see figure 14). ϵ is always about 2000.

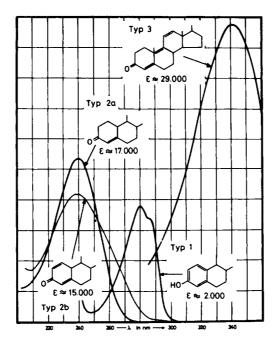


Fig. 14: Steroid hormones.

- 2. Among the 3-keto-steroids two types are distinguishable:
 - 2.1. the 4-en-3-one-steroids (progesterone [2208], testosterone [2220], norethisterone [2212], hydrocortisone [1008]) have a reasonably narrow absorption band at 240 nm, showing no absorption at wavelengths > 275 nm. ϵ is about 17000.
 - 2.2. the 1,4-diene-3-one steroids (prednisolone [1011], desoximetasone [1003], fluprednisolone [1006]) exhibit a maximum at 238 nm showing a broadening to the wavelength range up to 285. ϵ is about 15000 to 16000.

Otherwise the possibility of distinguishing various substances within these groups from their UV spectra is minimal and has to be restricted to ascertaining the E 1%/1cm (using pure substances).

Chemical, chromatographic and other spectroscopic methods must be employed for further identification.

3. Trien-one steroids (trenbolone 17 β -acetate [1014], norgestrienone [2207]) are unique: the long chromophore results in a longwave strong absorption band at 340 nm ($\epsilon \approx 29000$).

UV-Spektrokognosie

Es soll hier keine allgemeine Deutung von Spektren und keine Untersuchung aller denkbaren Struktureinflüsse gegeben werden – auch keine Beschreibung der Charakteristiken aller Stoffklassen. Ziel der nachstehenden Ausführungen, die auf UV-Spektren ausgerichtet sind, ist lediglich, dem nicht »hauptamtlichen« Spektroskopiker einige einfache Typen einzuprägen, die in der Arzneimittelanalyse häufig auftreten und die ihm dann in einfachen Fällen eine – wenigstens beschränkte – Zuordnung eines Spektrums zu einer bestimmten Stoffklasse ermöglichen soll. Daß komplexe Struktureinflüsse gelegentlich die sonst deutlichen und typischen Banden schwer erkenntlich oder unkenntlich machen können, muß bei der Deutung unbekannter Spektren immer berücksichtigt werden. Immerhin ist die Kenntnis wichtiger Grundtypen ein sehr wertvolles Rüstzeug für den Arzneimittelanalytiker.

Es soll gezeigt werden, wie man sich einige Grundkenntnisse zur UV-Spektrokognosie der Arzneistoffe aneignen und wie man die Kenntnis solcher »Familienmerkmale« verwerten kann. Manche Benutzer werden besondere Stoffklassen bearbeiten und ihre eigene Spektrokognosie entwickeln. Mit den folgenden Beispielen sollen in erster Linie Anregungen gegeben werden.

Phenyl - ohne sonstige chromophore Gruppen

Dieser Spektrentyp kommt im Bereich der Arzneistoffe häufig vor: Der Chromophor ist gekennzeichnet durch 3–5 relativ nahe zusammenliegende Banden im Bereich 246–268 nm. Die gute Auflösung der Feinbanden hängt von der Leistungsfähigkeit des Photometers, aber auch von der zu untersuchenden Substanz ab. Für eine gute Differenzierung sollte man darauf achten, daß in günstiger Konzentration gemessen wird. Die Wellenlängen der Banden sind in Abb.1 wiedergegeben.

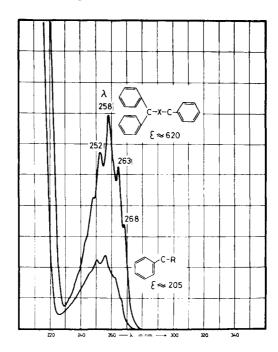


Abb. 1: Phenyl- ohne funktionelle Gruppen.

Die spezifische Absorption des reinen »Phenyl-Chromophors« ist recht niedrig: Für Verbindungen des Typs \bigcirc -R ist $\epsilon \approx 200$ (gemessen in Methanol). Treten 2 oder 3 Phenylreste in einem Molekül auf, so steigt ϵ auf etwa 400 bzw. 600. Entsprechend liegt die E 1%/1cm für viele Arzneistoffe etwa im Bereich 5–20, d. h., man sollte in Konzentrationen von ca. 20–50 mg/100 ml messen.

Zahlreiche Beispiele finden sich in der »Specfinder-Tabelle« im Bereich um 257 nm. Phenmetrazine hydrochloride [007] zeigt eine leichte hypsochrome Verschiebung nach 255 nm.

Beim benzyl penicillin [402] ist der kurzwellige Teil des Spektrums zwar »verwischt«, die Feinbanden bei 263 nm und 257 nm sind aber noch gut erkennbar.

Phenolische OH-Gruppen

	Spektrum	AbsMax.	Schulter	Veränderung
	Nr.	in Methanol	bei λ max.	in
			+8 nm	alk. Lösung
Phenol	1210	273 nm	+	7
Tyramine	2709	277 nm	+	7
Salicylic acid	1105	302 nm	Ø	Ľ
Salicylic acid methylester	703	305 nm	Ø	7
Salicylamide	118	288 nm	Ø	7
p-Hydroxybenzoic acid	1205	252 nm	Ø	7
p-Hydroxybenzoic acid ester	1207	256 nm	Ø	7
Paracetamol	112	247 nm	Ø	Z

Phenolische OH-Gruppen prägen das UV-Spektrum recht charakteristisch: Sind keine sonstigen funktionellen Gruppen am Benzolring vorhanden, so liegt das Absorptionsmaximum (in Methanol) bei etwa 274 \pm 3 nm mit einer sehr typischen Schulter bei 282 \pm 3 nm (vgl. Abb. 2). Hierher gehören neben phenol [1210], thymol [1212] und resorcinol [1104] viele Antihypotonika und Sympathomimetika der Gruppe 27. Bei den phenolischen Oestrogenen ist eine bathochrome Verschiebung um etwa 8 nm zu beachten (estradiol [2202], estradiol 17 β -cipionate [2203], estrone [2204], ethinylestradiol [2205]). Ortho-Diphenole (corbadrin [2701] usw.) zeigen die o. g. Schulter nicht. Ebenso nicht alle Verbindungen, die zusätzlich noch andere funktionelle Gruppen am Benzolkern haben. Bei solchen bi-funktionellen Verbindungen kommt es stets zu starken Verschiebungen des gesamten Spektrums zum Lang- oder Kurzwelligen hin.

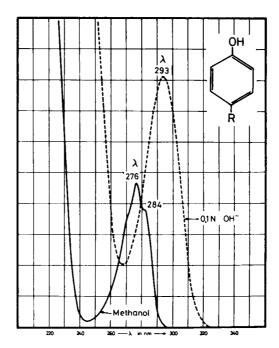


Abb. 2: Phenol. OH- ohne sonstige funktionelle Gruppen.

Fast alle phenolischen Verbindungen zeigen in wäßrig-alkalischer Lösung – gegenüber Methanol – infolge Dissoziation eine sehr deutliche bathochrome Verschiebung, meist (aber nicht immer) verbunden mit einer Erhöhung der Absorption (Symbol: 7). Salicylic acid [1105] macht eine Ausnahme: Hier überwiegt der Einfluß der freien ortho-ständigen Carboxylgruppe. Bei salicylic acid esters [703] tritt die Alkali-Verschiebung wieder auf. Ortho-Diphenole sind wegen ihrer Zersetzlichkeit in alkalischer Lösung kaum meßbar.

Achtung: Bathochrome Verschiebung → oder ७ in alkalischer Lösung tritt nicht ausschließlich bei Phenolen auf! Sie ist auch bei hydroxypyridines [2806] und enolisierbaren Ketonen zu beobachten: pyridone [909], benzoxazolinone [2402], bei benzodiazepinones [2604], außerdem bei Nitroverbindungen (metronidazole [916], nifurprazine [917], nitrofural [918], nitrofurantoin [919]).

UV-Spektrokognosie Fortsetzung

Immerhin kann man die Schulter bei λ max. + 8 nm zusammen mit der π -Verschiebung als guten Hinweis auf phenolisches OH nehmen – insbesondere bei Abwesenheit anderer funktioneller Gruppen am Ring. Die π -Verschiebung allein gilt nur als Hinweis mit Einschränkungen.

Die Intensität der Phenolbande bei ca. 275 nm ist mittelstark: ϵ liegt in der Größenordnung um 2000 (abhängig auch von sonstigen Substituenten), bei Ortho-Diphenolen um 3000.

Phenolether

Phenolether, insbesondere solche ohne weitere funktionelle Gruppen, zeigen in methanolischer Lösung ein relativ ähnliches UV-Spektrum wie Phenole: λ max. etwa 270–285 nm, teils ohne deutliche Schulter (meticillin sodium salt [408], verapamil hydrochloride [1416], teils mit schwacher Schulter (guaifenesin [2406]), zum Teil mit deutlicherer Feinbande bei λ max. + 8 nm (chlorphenesin carbamate [1101]).

Eine bathochrome Verschiebung in alkalischer Lösung findet hier natürlich nicht statt (Unterscheidung morphine hydrochloride [111] – codeine phosphate [801]). Bei vielen anderen Verbindungen ist die nur relativ schwache Phenoletherbande stark überlagert und kaum mehr zu erkennen (mepyramine maleate [306], noscapine hydrochloride [805], papaverine hydrochloride [2408]).

Die Intensität der Phenolether-Absorption ist niedriger als die der Phenole: ϵ liegt in der Größenordnung um 1000.

Aromatische Carbonsäuren

Aromatische Carbonsäuren (benzoic acid [2902], aspirin [102]) stellen einen mäßig charakteristischen Chromophor dar: Absorptionsmaxima liegen bei etwa 273 \pm 3 nm und etwa 227 \pm 3 nm und sind nur zusammen als Hinweis zu bewerten. Durch Dissoziation in NaOH wird die Absorption gegenüber der Lösung in Methanol geschwächt. Gleichzeitig vorhandene andere Chromophore schlagen meist durch.

Bei den Pyridincarbonsäuren tritt nur bei der isonicotinic acid [2907] die Carboxylfunktion im Spektrum deutlicher mit in Erscheinung.

Aromatische Amine

Aromatische Amine zählen zu den stärksten und für die Arzneimittelanalytik wichtigsten Chromophoren. Ihr Absorptionsmaximum liegt um 290 nm, wird aber durch verschiedene Substituenten und funktionelle Gruppen zum Teil auch erheblich verschoben. Die hohe Intensität (in Methanol) ist durch $\epsilon \approx 20000$ gekennzeichnet; daher schlägt dieser Chromophor praktisch immer sehr deutlich durch. Außerordentlich charakteristisch ist die Unterdrückung dieser Bande – praktisch bis auf 0 – bei zunehmender H $^+$ -Konzentration der Lösung, d. h. zunehmender Dissoziation (vgl. Abb. 3). Man kann auf dieser Basis den pKa-Wert aromatischer Amine bestimmen. Soweit im kurzwelligen Bereich keine Überdeckung durch andere Molekülkomponenten stattfindet, bildet sich mit zunehmender Dissoziation eine Bande um etwa 220–240 nm neu oder verstärkt aus (benzocaine [1701], tetracaine hydrochloride [1709], acriflavinium chloride [902], sulfaguanidine [928]). In 0,1 M Natriumhydroxid wird die Hauptbande im allgemeinen leicht zum Kurzwelligen hin verschoben.

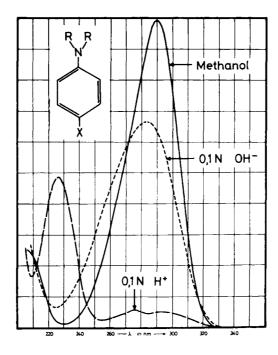


Abb. 3: Aromatische Amine.

Primäre aromatische Amine zeigen im Prinzip gleichartige Absorptionsmerkmale wie sekundäre und tertiäre (mit aliphatischen Resten).

Sulfonamide

Sulfonamide (phthalylsulfathiazole [920], sulfacetamide [921], sulfaquinoxaline [922] usw.) mit chemotherapeutischer Wirkung sind in allen Fällen 4-Amino-benzolsulfonamide (Sulfanilamide). Dementsprechend ist ihre UV-Absorption vorwiegend durch die aromatische Aminogruppe bestimmt. Sulfanilamide zeigen eine Hauptbande meist bei etwa 270 nm.

Sie liegt also deutlich kurzwelliger als z. B. bei den 4-Aminobenzoesäureestern. Die Abschwächung der Absorption in saurer Lösung ist nicht überall gleich stark ausgeprägt.

Selbstverständlich wird die Lage der Hauptbande und auch die Veränderung infolge Dissoziation durch die Substituenten sowohl an der SO_2 -NH-Gruppe als auch an der Aminogruppe beeinflußt: Starke chromophore Komponenten wie Phthalyl- oder Chinoxalyl- [920, 922] u. a. schlagen dementsprechend deutlich durch. Aus diesem Grund sind die UV-Spektren chemotherapeutischer Sulfonamide nicht in allen Fällen auf den ersten Blick als »typisch« zu erkennen.

Sulfonamid-Diuretika

Sulfonamid-Diuretika (cyclothiazide [1307], furosemide [1309], hydrochlorothiazide [1310], bendroflumethiazide [1315], hydroflumethiazide [1316], butizide [1317], cyclopenthiazide [1318]) der allgemeinen Formel

erkennt man leicht an einer starken Absorptionsbande bei 270–278 nm ($\epsilon \approx 20000$), einer schwächeren Bande bei 315–325 nm ($\epsilon \approx 3000$ –4000) und einem tiefen Minimum bei 240 nm (vgl. Abb. 4). Bei furosemide [1309] ist die langwellige Bande unter dem Einfluß des Furylsubstituenten weiter bathochrom (\rightarrow 336 nm) verschoben.

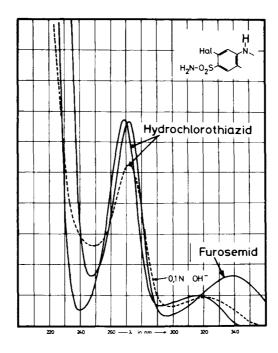


Abb. 4: Sulfonamid-Diuretika.

Chlorothiazide [1305] (ohne N-**H** am Benzolsulfonamidring), bumetanide [1302] und chlortalidone [1303] gehören nicht zum o. g. Spektrentyp.

Anilide

Säureamide aromatischer Amine müssen für die Arzneimittel-Spektrokognosie in zwei weitgehend verschiedene Gruppen unterteilt werden:

- A) die der einfachen Säureanilide vom Typ acetanilide/phenacetine [101, 114] und
- B) die am aromatischen Ring in 2,6-Stellung substituierten Anilide vom Typ lidocaine/butanilicaine [1707, 1703]

Typ A zeigt eine starke Absorptionsbande bei 240–250 nm (ϵ etwa 1500), die in Methanol deutlich stärker ist als in wäßriger Lauge oder Säure (vgl. Abb. 5). Auch paracetamol [112] gehört in diese Gruppe; es läßt sich von phenacetin [114] durch eine $\mbox{1}$ -Verschiebung in 0,1 M OH $^{-}$ (infolge Dissoziation der Phenolgruppe) gut unterscheiden.

Einigermaßen charakteristisch für phenacetin [114], bucetin [104] und paracetamol [112] ist eine leichte Schulter bei 280–290 nm.

Anmerkung: Acetohexamide [501] zeigt – trotz anderer Struktur – ein sehr ähnliches Spektrum.

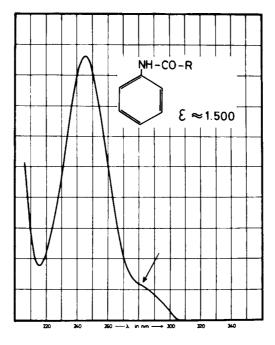


Abb. 5: Acetanilide Typ A.

Die 2,6-substituierten Anilide Typ B zeigen nur eine schwache Bande um 265 nm ($\epsilon \approx 350$) mit einer Feinbande bei λ max. + 8 nm (vgl. Abb. 6).

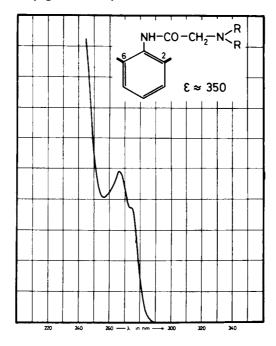


Abb. 6: Acetanilide Typ B.

In wäßrigen Lösungen ist hier die Absorption gegenüber der Lösung in Methanol leicht erhöht (↑). Im Bereich um 245 nm zeigen die Verbindungen vom Typ B nur eine unspezifisch zum kurzwelligen UV hin ansteigende Absorption, aber keine ausgeprägte Bande wie Typ A.

Pyrazolone

Die Pyrazolone des Arzneimittelbereichs stellen UV-spektroskopisch eine relativ geschlossene und leicht erkennbare Gruppe dar, zu der auch das phenylbutazone [701] hinzugenommen werden kann. Typisch sind 2 relativ breite (unterschiedlich stark ausgeprägte) Maxima, die sich in ihrer Lage und Höhe etwas unterscheiden: Sie liegen um 234–246 nm und 265–275 nm (vgl. Abb. 7). Bei den basisch substituierten Verbindungen (dimethylamino phenazone [103], metamizole sodium [109]) tritt in saurer Lösung das langwellige Maximum etwa an die Stelle des Minimums in methanolischer Lösung. In alkalischer Lösung sind die Kurven etwas abgeflacht. ϵ liegt in methanolischer Lösung im Bereich 10000.

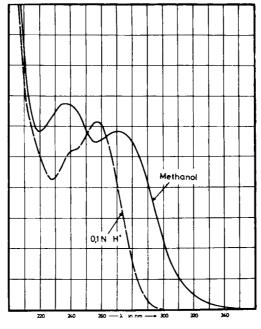


Abb. 7: Pyrazolone.

Pyridyl

Die Absorption des Pyridylrests ist wesentlich stärker als die des unsubstituierten Phenyl-. Er zeigt eine scharfe Bande bei etwa 261 nm mit Kleinbanden oder Schultern bei 269 und 265 nm (vgl. Abb. 8).

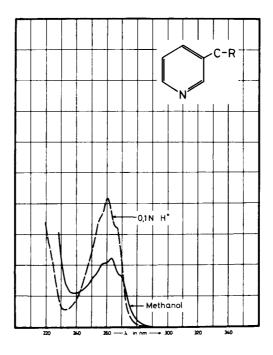


Abb. 8: Pyridyl-.

Ein wichtiges Merkmal ist die starke Erhöhung der Bande in 0,1 M Säure (Beispiele: pheniramine maleate [308], nicotine [3003]). Diese Charakteristika bleiben zum Teil auch in Gegenwart funktioneller Gruppen erhalten, so bei nicotinic acid und nicotinic acid amide [1802, 2805]. Bei isonicotinic acid [2907] gehen die Feinbanden verloren; bei hydroxy-pyridine-Verbindungen [2806] verändert sich der Charakter des Spektrums völlig, eingeschränkt vergleichbar mit den Phenolen (z. B. →-Verschiebung in alkalischer Lösung).

Indole

Indole zeigen ein recht charakteristisches, leicht zu erkennendes Spektrum mit einem Gipfel bei 280 nm, einer Schulter bei ca. 274 nm und einer scharfen Feinbande bei 289 \pm 1 nm (vgl. Abb. 9). »Schulbeispiele« aus dieser Sammlung sind tryptamine [2909] und tryptophan [3006].

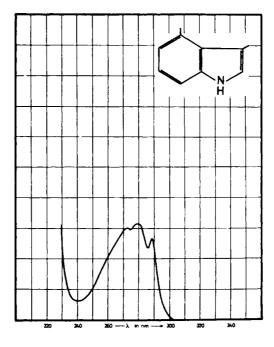


Abb. 9: Indolyl.

Diese Indol-Charakteristika treten auch bei den hydrierten Mutterkornalkaloiden, nicht mehr aber bei den Lysergsäurederivaten in Erscheinung.

Bei pindolol [2008] ist die Feinbande bei 288 nm noch sehr deutlich; der Indol-Gipfel bei 280 nm tritt nur noch als Schulter in Erscheinung (das Maximum von pindolol bei 265 nm ist von der Indolether-Gruppe mit geprägt).

Aus solchen Beispielen wird deutlich, daß man sich beim Bestreben um Zuordnung nicht nur das Maximum ansehen sollte.

Xanthine

Die Xanthine des Arzneischatzes stellen sich UV-spektrophotometrisch als eine sehr homogene Gruppe dar: Alle hier aufgeführten Verbindungen (caffeine [002], pentoxifylline [006], diprophylline [1406], etofylline [1407], theobromine [1414], theophylline [1415]) zeigen eine schöne und konstante Absorptionsbande bei 272 \pm 2 nm mit einem ϵ von ca. 9000–10000 (in Methanol): Minimum bei 242 \pm 2 nm (vgl. Abb. 10).

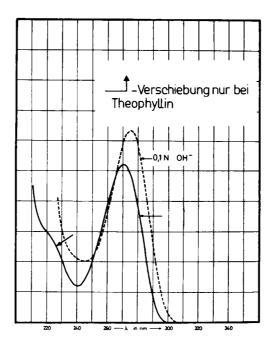


Abb. 10: Xanthine.

Die Erkennung der Xanthine wird gut abgesichert durch eine schwache, aber sehr typische Schulter bei 225–230 nm. Nur theophylline [1415], die einzige nennenswert dissoziierbare Verbindung dieser Reihe, gibt eine einigermaßen deutliche ¬Verschiebung in Alkalilauge.

Phenothiazine

Phenothiazine sind aus dem UV-Spektrum leicht als solche zu erkennen: Eine breitere, niedrigere Bande im Langweiligen und eine scharfe Bande bei etwa 250 nm (ϵ ca. 30000) kennzeichnen den Verbindungstyp zusammen mit den Minima (vgl. Abb. 11) .

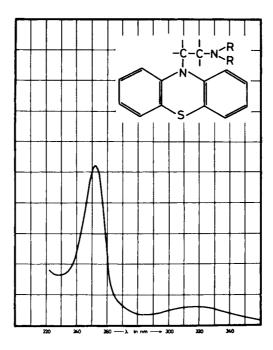


Abb. 11: Phenothiazine.

Die Neuroleptika mit basischer Seitenkette (chlorpromazine hydrochloride [2102], profenamine hydrochloride [2107]) zeigen gegenüber dem Grundkörper (phenothiazine [203]) eine leichte Verschiebung zum Kurzwelligen hin.

Benzodiazepine

Benzodiazepine der bekannten Reihe 1,4-Benzodiazepin-2-on zeigen in saurer Lösung eine charakteristische, leicht asymmetrische Absorptionsbande bei etwa 283 nm ($\epsilon \approx 12000$) sowie ein flacheres Maximum bei 360 nm (vgl. Abb. 12). Diazepam [2604] läßt sich von oxazepam [2605] dadurch unterscheiden, daß bei letzterem eine deutliche Rotverschiebung der flachen Bande auftritt (in Methanol 315, in Alkali \rightarrow 344 nm).

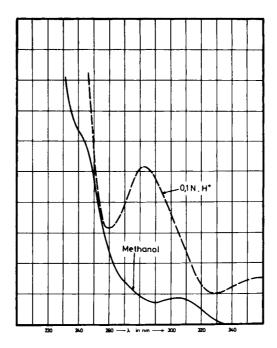


Abb. 12: 1,4-Benzodiazepin-2-one.

Die Spektren von chlordiazepoxide hydrochloride [2602] und von clobazam [2603] (1,5-Benzodiazepin-2,4-dion) sind von denjenigen der vorgenannten Reihe deutlich unterschieden.

Barbiturate

Barbiturate sind UV-spektrophotometrisch wenig ergiebig. Sie zeigen nur in alkalischer Lösung einigermaßen ausgeprägte UV-Spektren im ziemlich kurzwelligen Bereich. Gegenüber den rein aliphatisch substituierten Vertretern (λ max. um 242 nm) zeigt phenobarbital [1508] ein Maximum bei 252 nm infolge des verlängerten Chromophors (vgl. Abb. 13).

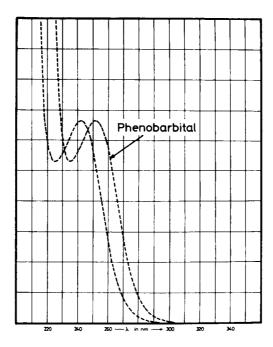


Abb. 13: Barbiturate in 0,1 M OH⁻.

Steroidhormone

Steroidhormone lassen sich nach ihren UV-Spektren leicht in 3 Gruppen einteilen:

1. Die Oestrogene (estradiol 17β -cipionate [2203], estrone [2204], ethinylestradiol [2205]) mit ihrer phenolischen OH-Gruppe zeigen ein Absorptionsverhalten, wie zuvor für Phenole beschrieben; die Spektren sind gegenüber einfachen Phenolen lediglich leicht in den längerwelligen Bereich verschoben (λ max. um 280 nm) (vgl. Abb. 14). ϵ ist jeweils etwa 2000.

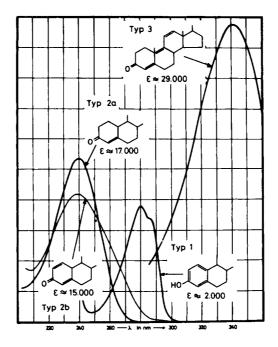


Abb. 14: Steroidhormone.

- 2. Bei den 3-Keto-Steroiden lassen sich bei genauer Betrachtung (!) noch 2 Typen unterscheiden:
 - 2.1. die 4-en-3-on-Steroide (progesterone [2208], testosterone [2220], norethisterone [2212], hydrocortisone [1008]) zeigen eine (relativ) schmale Bande mit Maximum bei 240 nm, die oberhalb 275 nm keine Absorption mehr aufweist; $\epsilon \approx 17000$.
 - 2.2. die 1,4-dien-3-on-Steroide (prednisolone [1011], desoximetasone [1003], fluprednisolone [1006]) zeigen ein Maximum bei 238 nm mit einem zum Langwelligen (bis 275 nm) hin verbreiterten Bandenfuß; $\epsilon \approx 15000-16000$.

Im übrigen ist das Differenzierungsvermögen durch UV innerhalb dieser Gruppen minimal und muß sich (bei Einsatz reiner Substanzen!) auf eine Ermittlung der E 1%/1cm beschränken. Hier müssen chemische, chromatographische und andere spektroskopische Methoden weiterhelfen.

3. Eine Sonderstellung nehmen die (Ges)-trienone (trenbolone 17β -acetate [1014], norgestrienone [2207]) ein: Der hier vorliegende lange Chromophor erbringt eine langwellige starke Absorptionsbande bei 340 nm ($\epsilon \approx 29000$).

Specfinder


- In the Specfinder table you can search for UV absorption maximum values. This is possible by entering the particular figure + nm (without blank; e. g. 348nm) in the search function of Acrobat Reader. The found spectrum can be called up by clicking on the respective spectrum number.
- In der Specfinder-Tabelle können Sie nach UV-Absorptionsmaxima suchen. Geben Sie dazu in der Suchfunktion des Acrobat Reader die gewünschte Zahl + Angabe nm ein (ohne Leerschritt; z.B.: 348nm). Das gefundene Spektrum läßt sich durch Anklicken der entsprechenden Spektren-Nummer aufrufen.

© 2002 Editio Cantor Verlag Aulendorf

Further tables

- Pharmaceutical substances showing a well distinct UV spectrum in alkaline solution or being insoluble in methanol
- Pharmaceutical substances showing no distinct <u>maximum</u> of absorption
- Pharmaceutical substances without useful UV absorption >220 nm



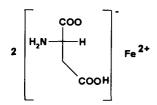



 M_{r}

320.0


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	Decom- position observed		
E ^{1%} 1cm				
ε				

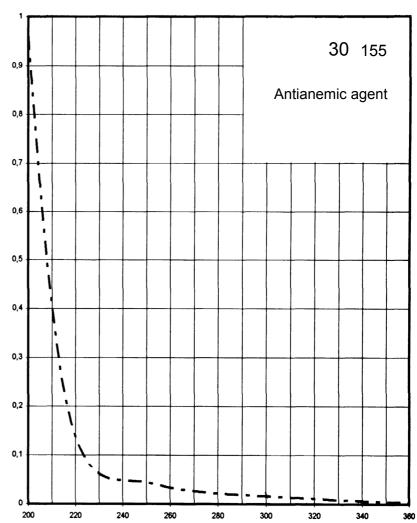

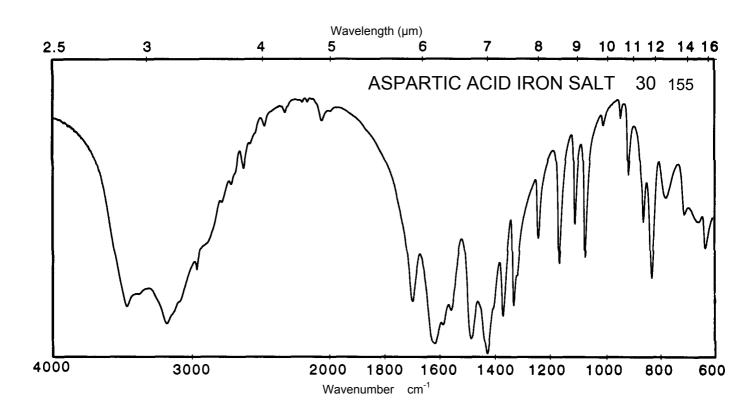



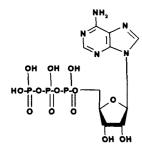



 M_{r}

507.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm	259 nm	257 nm	260 nm
E ^{1%} 1cm	247	247	239	251
3	12500	12500	12100	12700

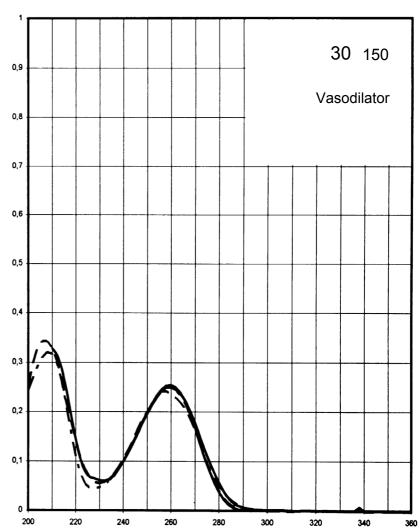

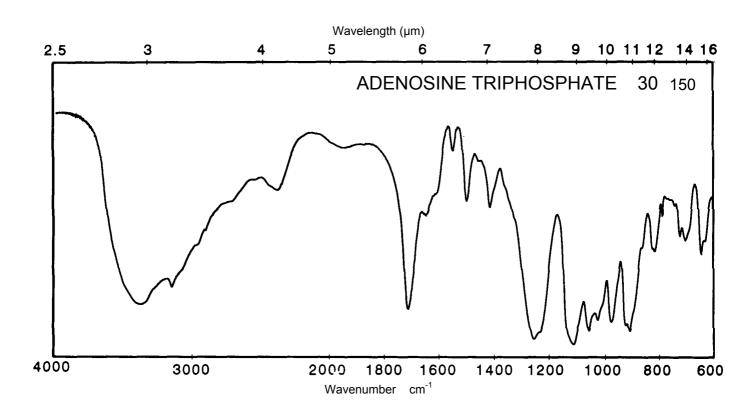

Name ADENOSINE PHOSPHATE

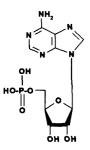



 M_{r}

347.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm	259 nm	257 nm	260 nm
E ^{1%} 1cm	414	397	398	401
3	14400	13800	13800	13900

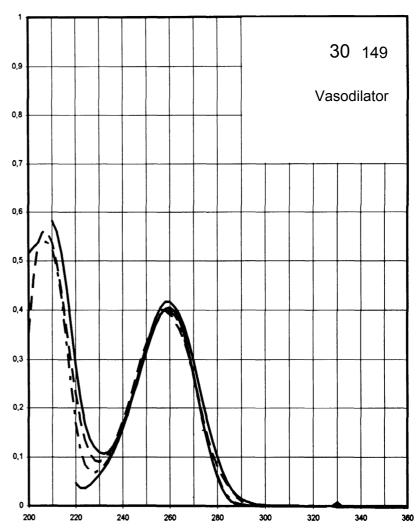

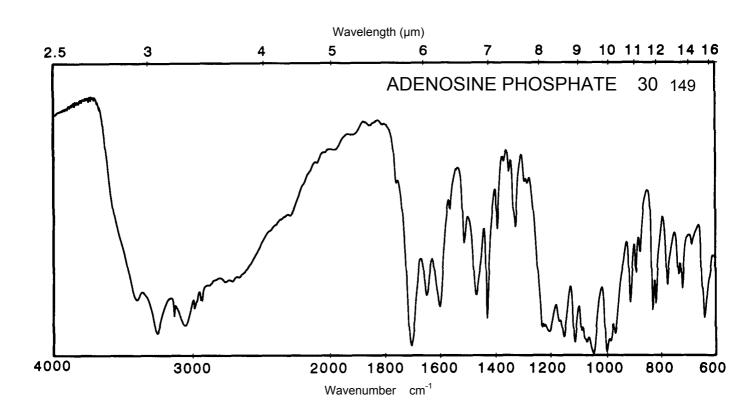

Name ACTINOQUINOL SODIUM SALT

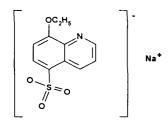


 M_{r}


275.3

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	310 nm 240 nm	253 nm	305 nm 240 nm
E ^{1%} 1cm		190 1200	1365	196 1205
3		5200 33100	37600	5400 33100

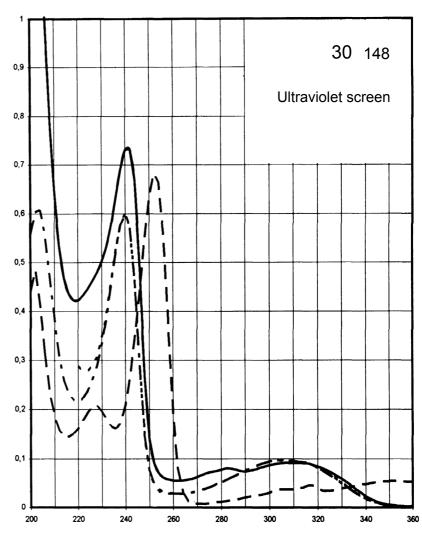

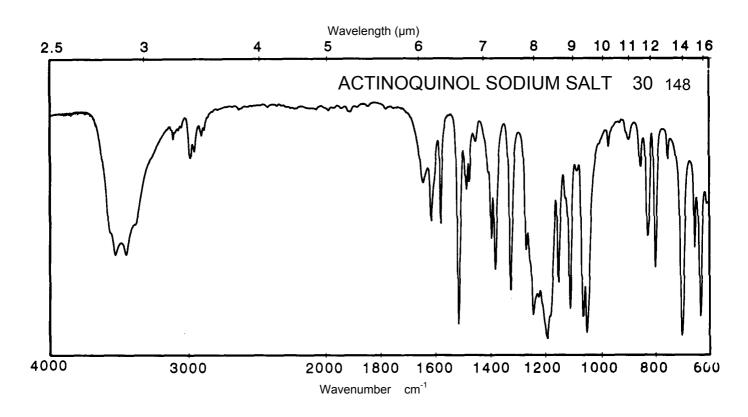



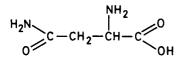



 M_{r}

132.1


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

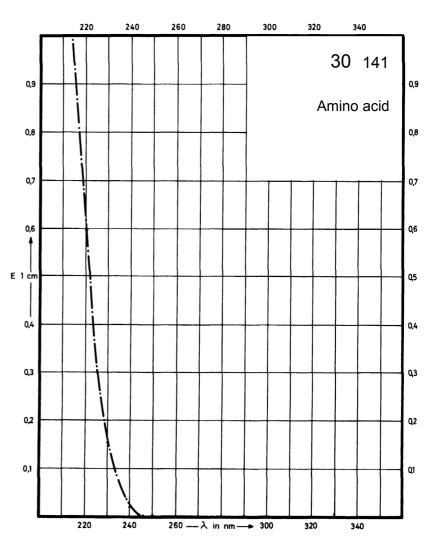

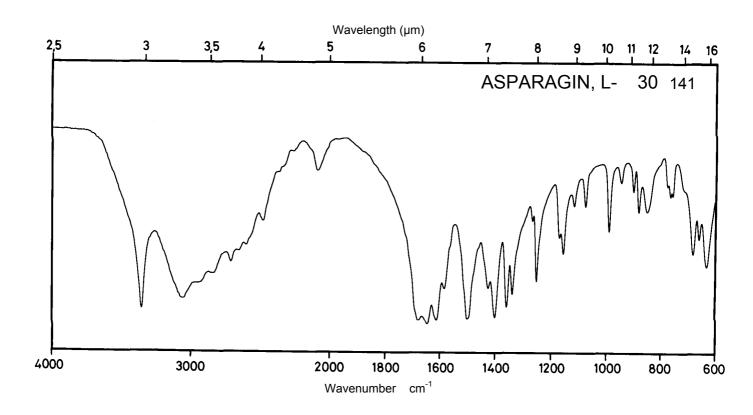



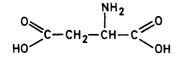



 M_r

133.1


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

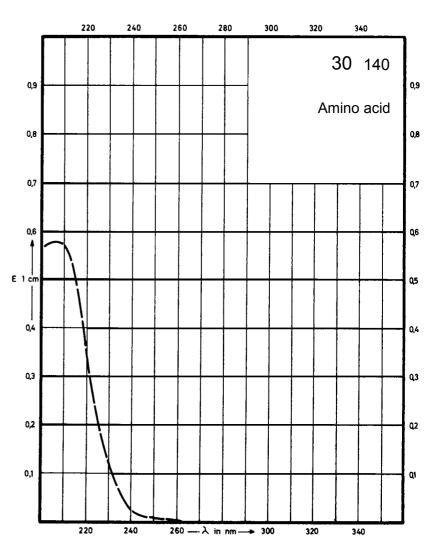

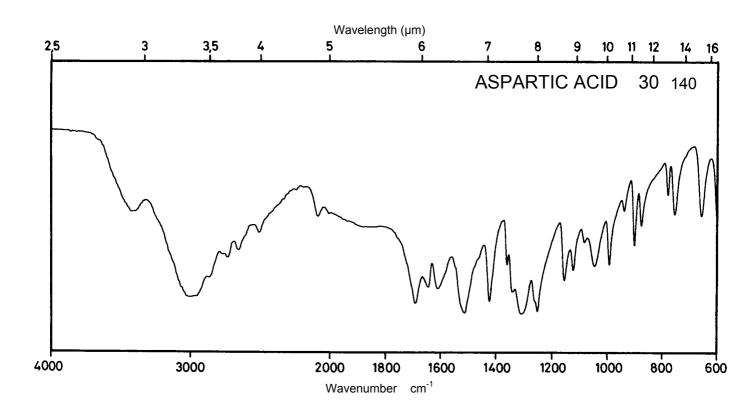




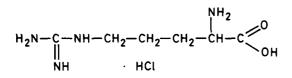



 M_r 210.7


Concentration 50 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

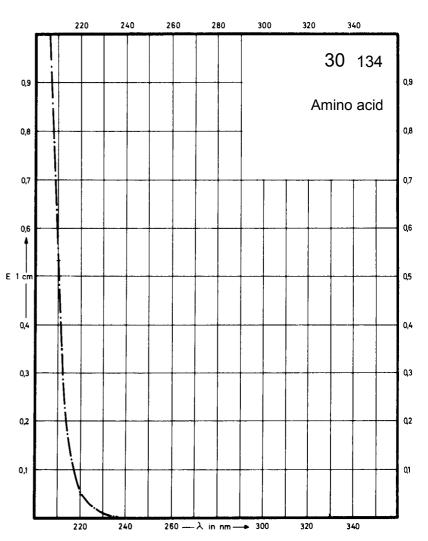


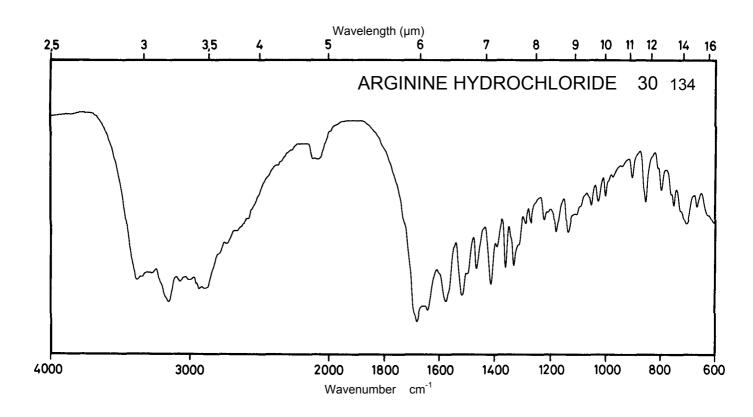

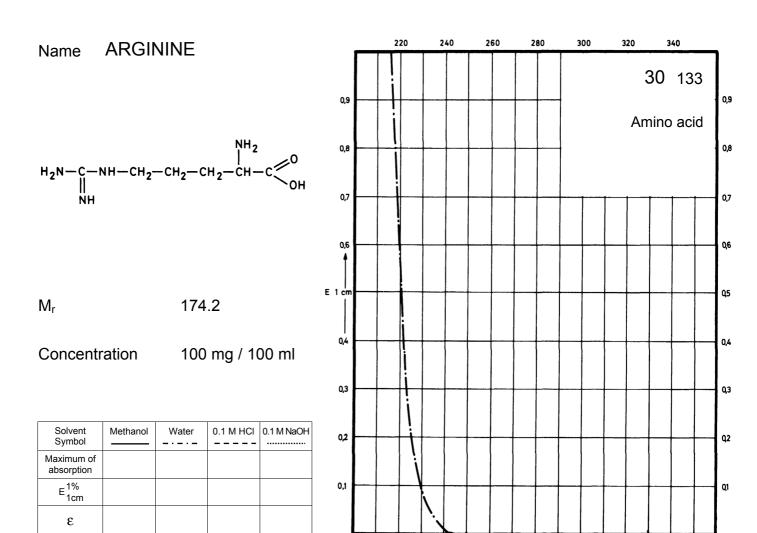


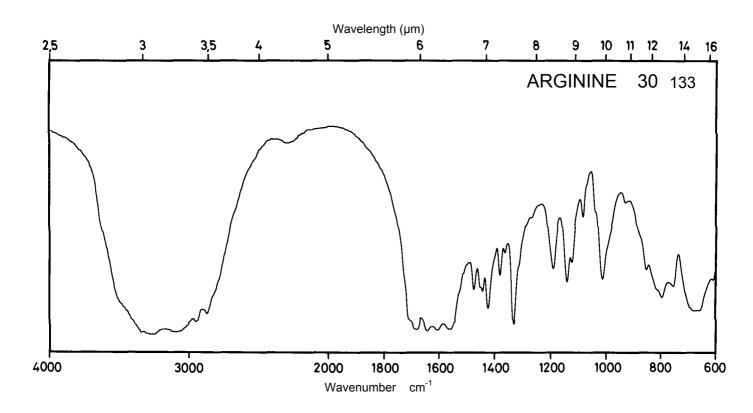

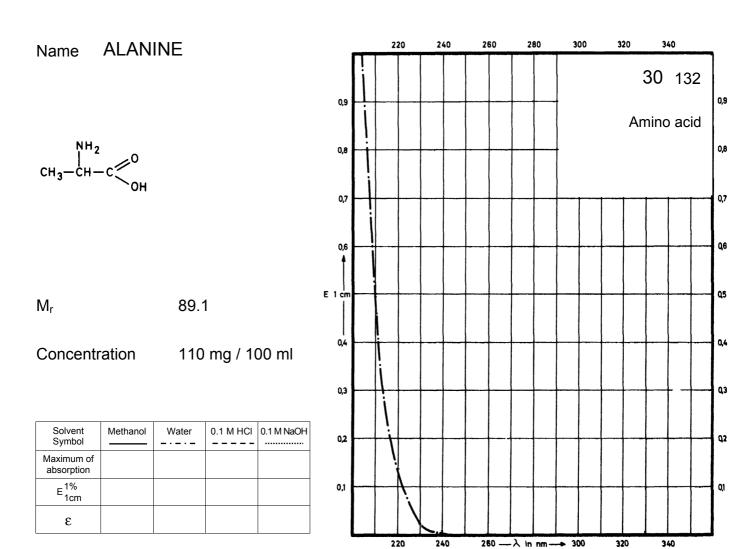

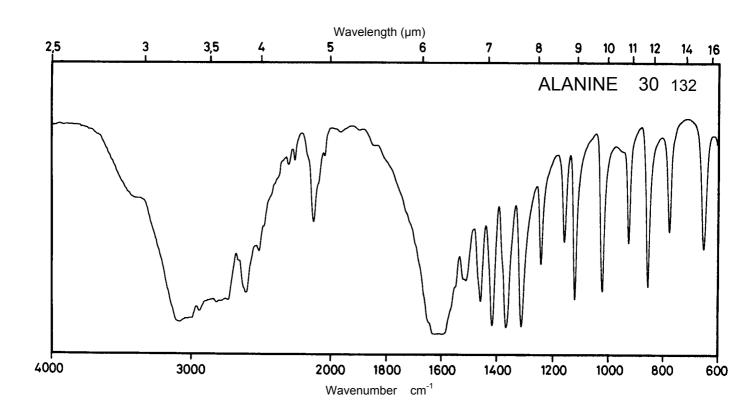



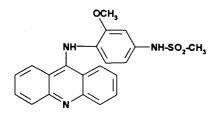



393.5


Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	434 nm 264 nm 248 nm	435 nm 264 nm	Decom- position observed	Decom- position observed
E ^{1%} 1cm	288 1020 1160	330 1326		
ε	11300 40200 45500	13000 52200		

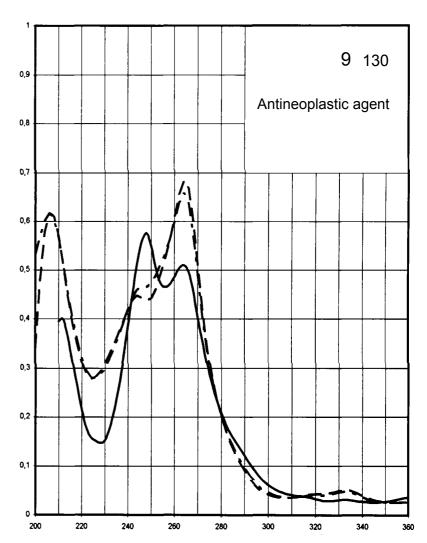


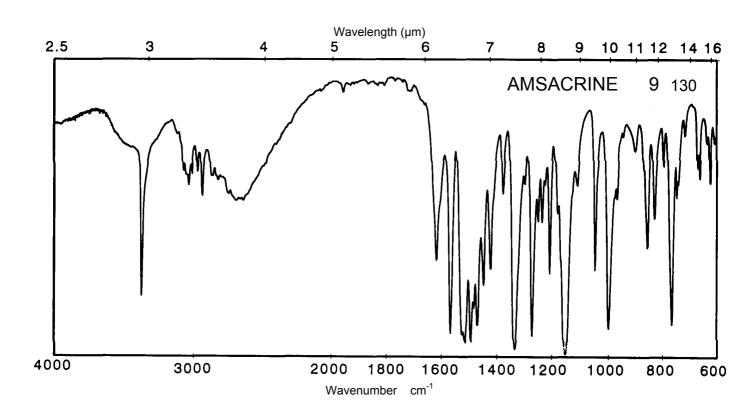

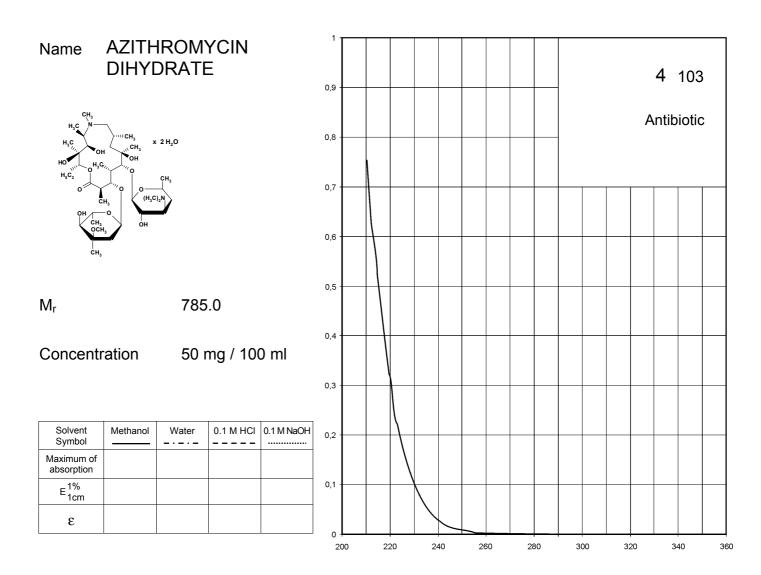

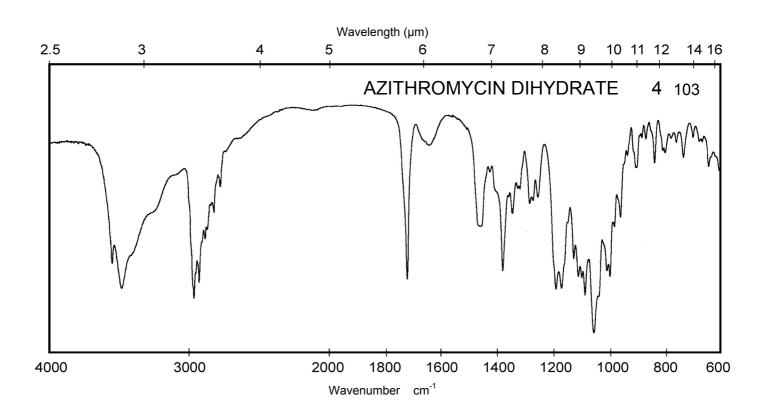





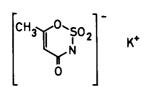


 M_r 201.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	228 nm	226 nm	225 nm	
E ^{1%} 1cm	543	557	493	
ε	10920	11200	9920	

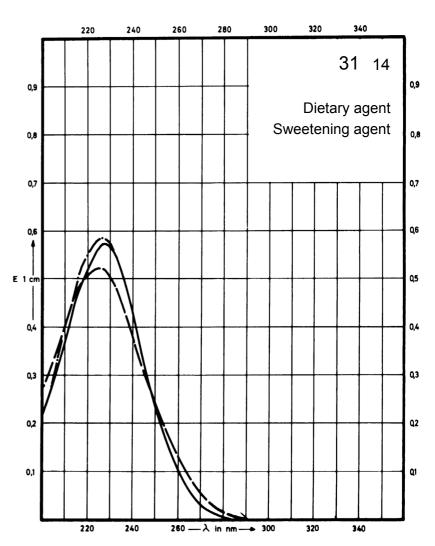

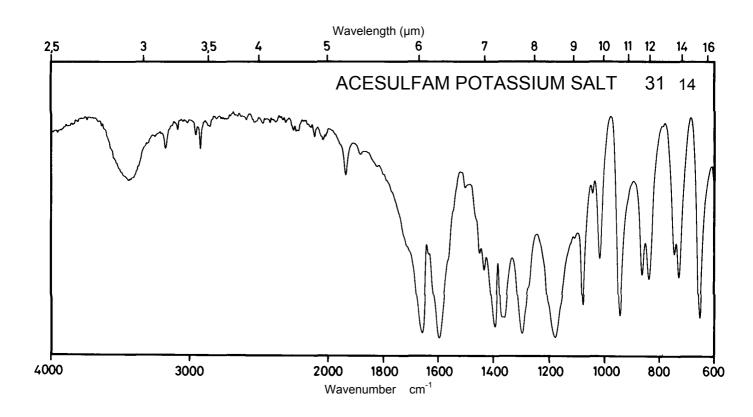





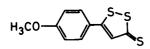



240.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	428 nm 348 nm 233 nm			
E ^{1%} 1cm	491 858 552			
ε	11800 20600 13300			

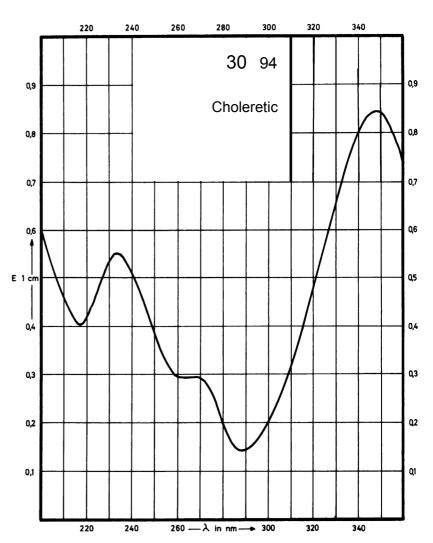


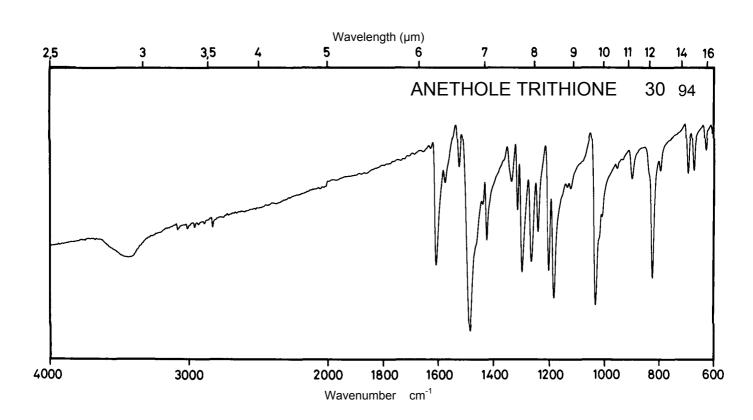

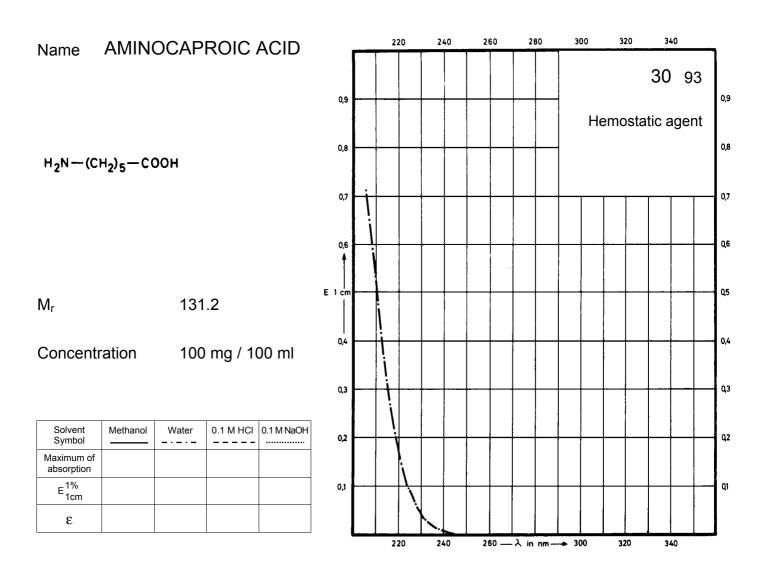

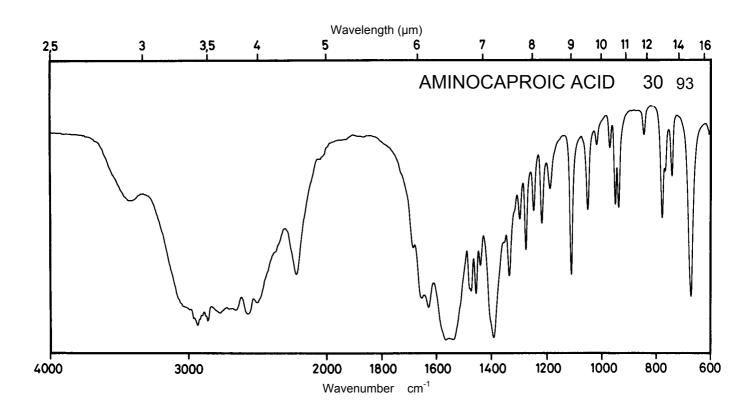





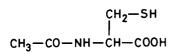


163.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

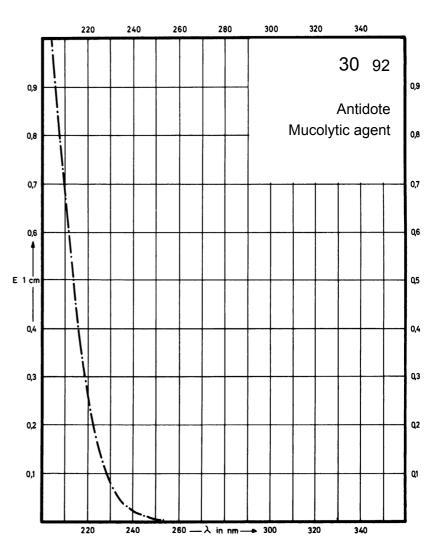

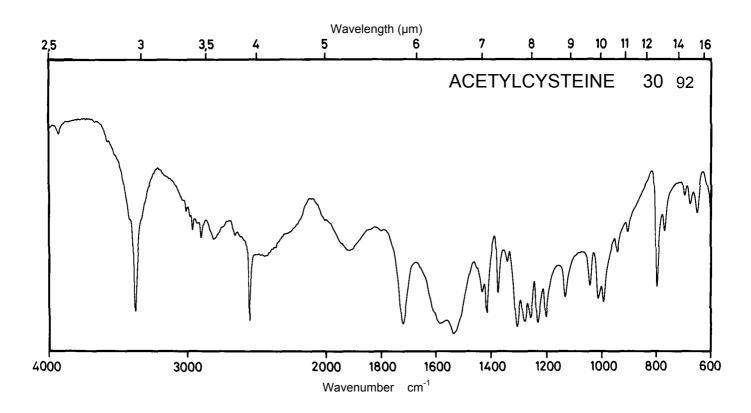



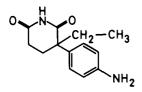



232.3

Concentration 1


1.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	242 nm			233 nm
E ^{1%} 1cm	500			Decom- position observed
3	11600			

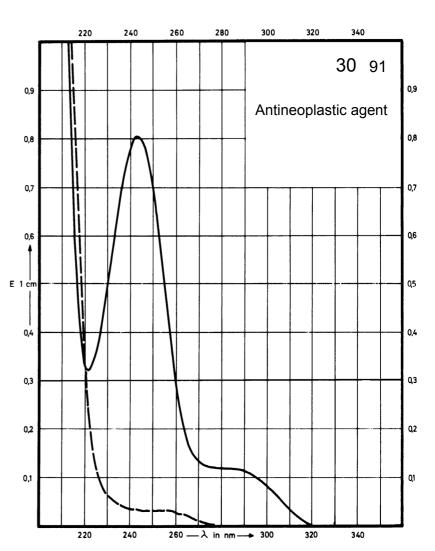

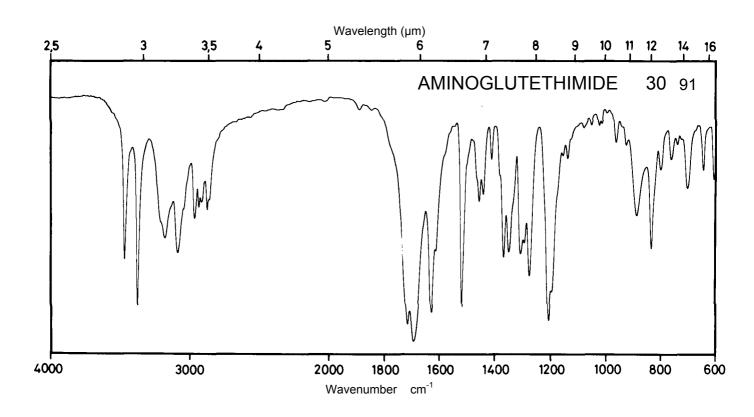





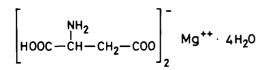


 M_{r} 360.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

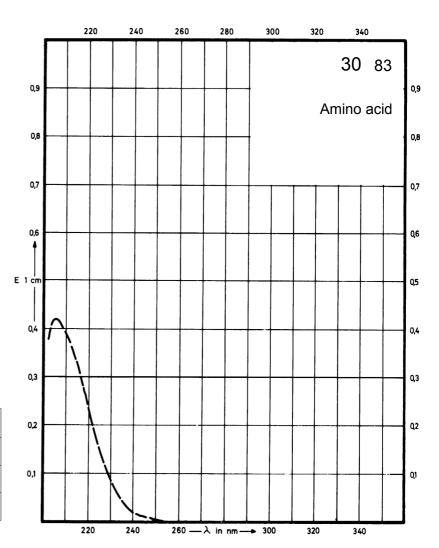

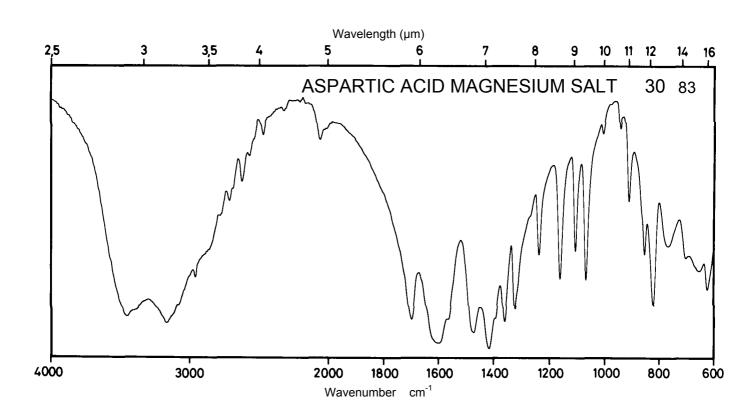




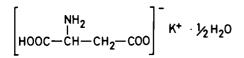


 M_r

180.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

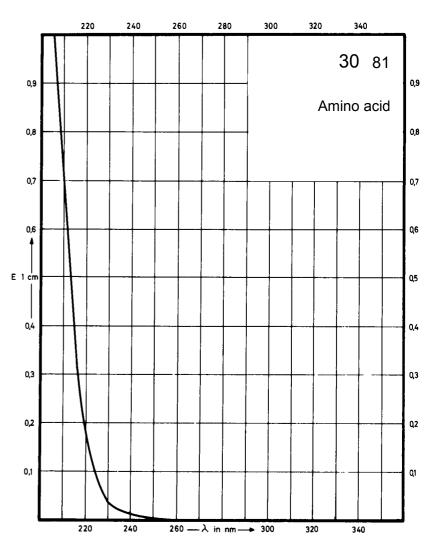

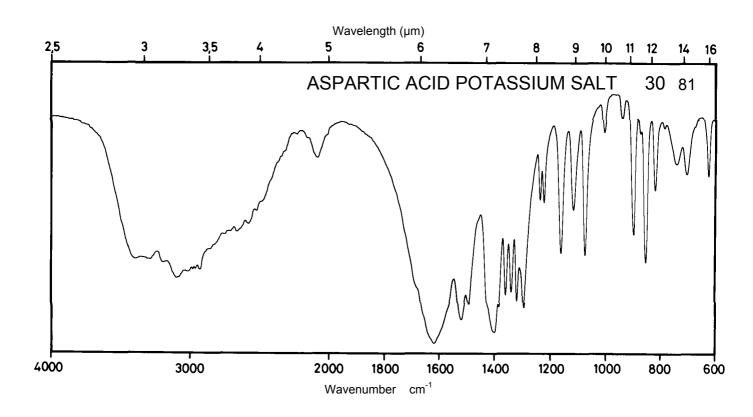



Name AZATHIOPRINE


H₃C - N NO₂


 M_r 277.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm		280 nm	285 nm
E ^{1%} 1cm	633		621	560
ε	17600		17200	15500

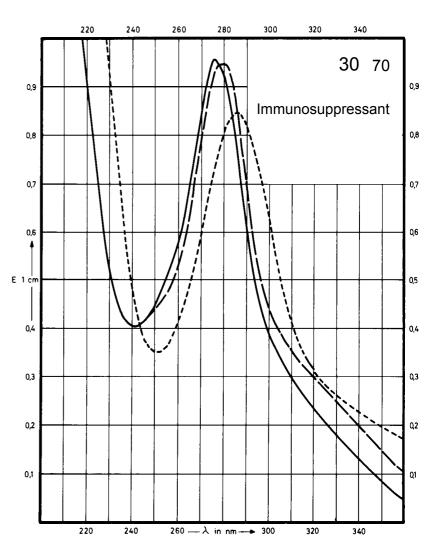

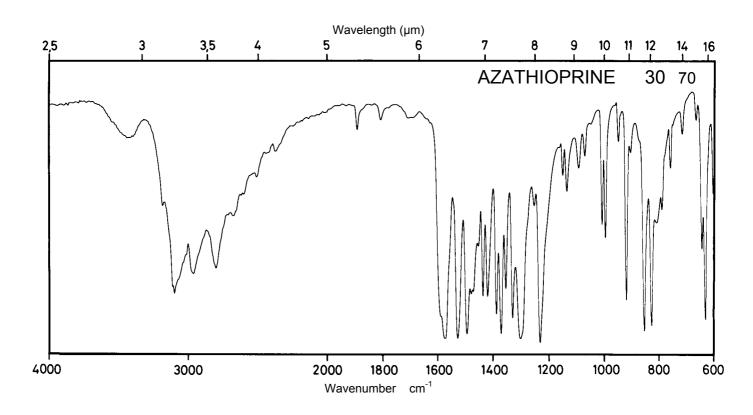

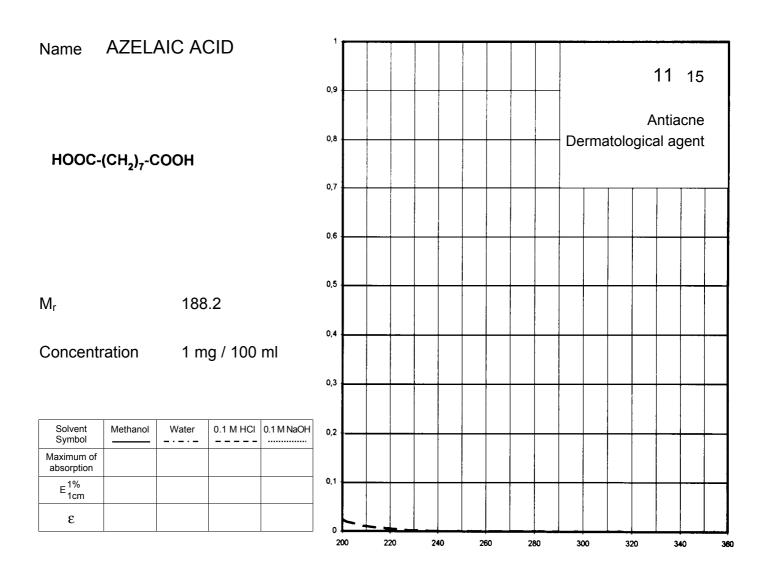

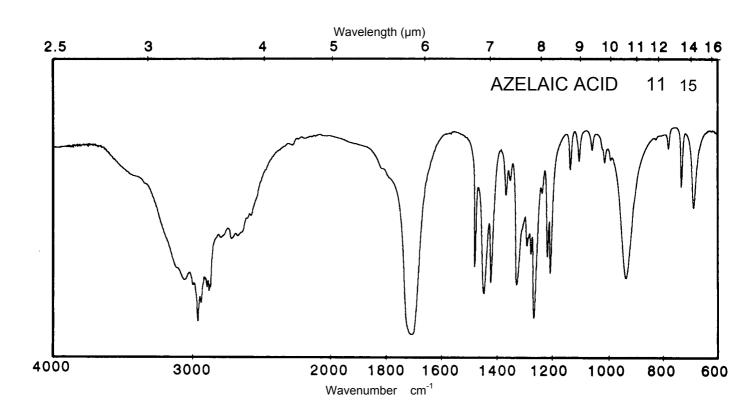



Name AMFEPRAMONE HYDROCHLORIDE

 M_r


241.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	249 nm	252 nm	252 nm	246 nm
E ^{1%} 1cm	563	565	560	462
3	13620	13660	13540	11170

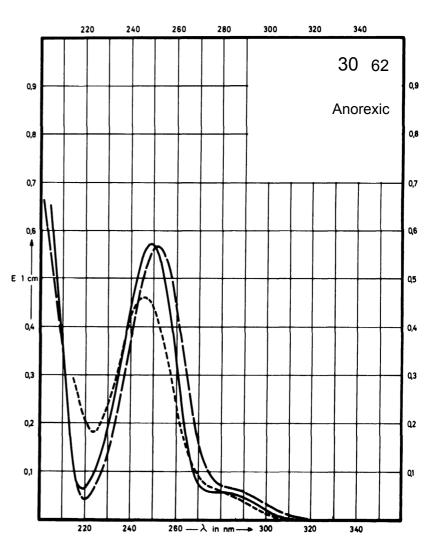

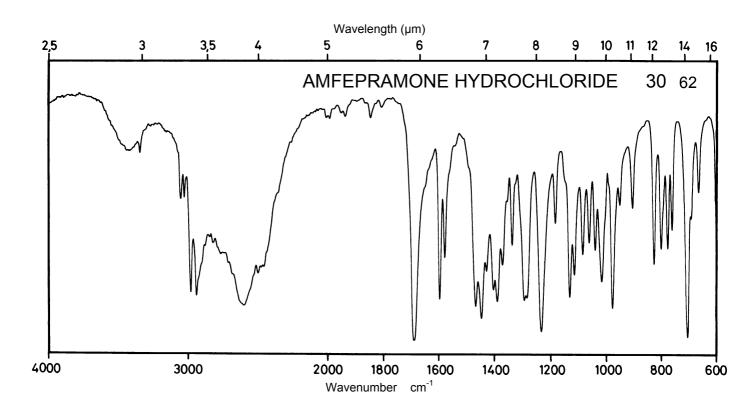





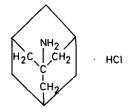



187.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

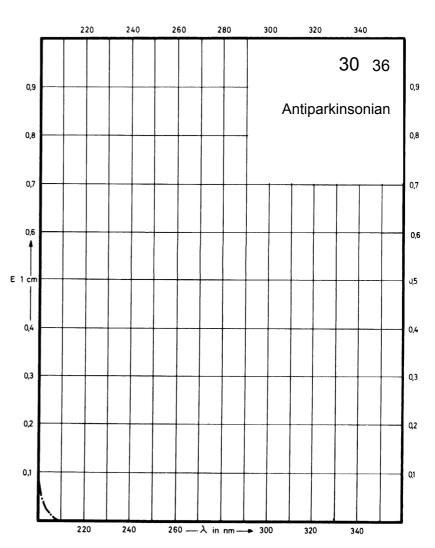

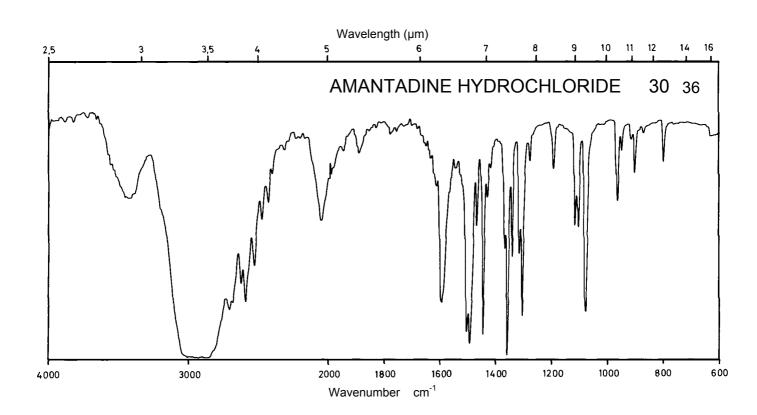


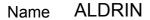


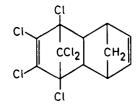

 M_r 364.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

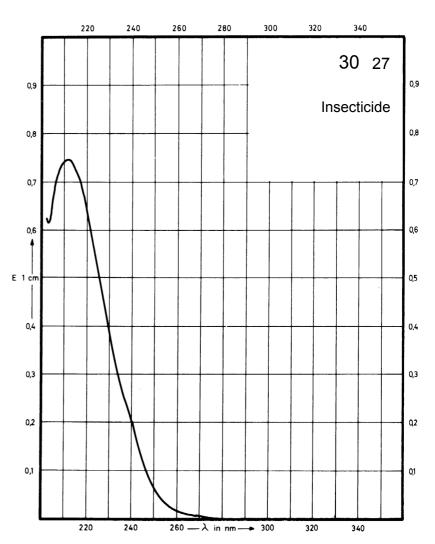

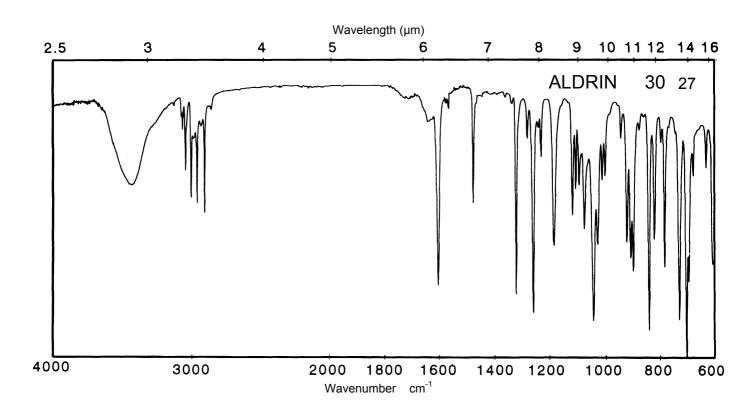



Name ACENOCOUMAROL

 M_{r}


353.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	305 nm 282 nm		304 nm 284 nm	300 nm
E 1% 1cm	510 637		573 637	598
ε	18040 22510		20240 22510	21120

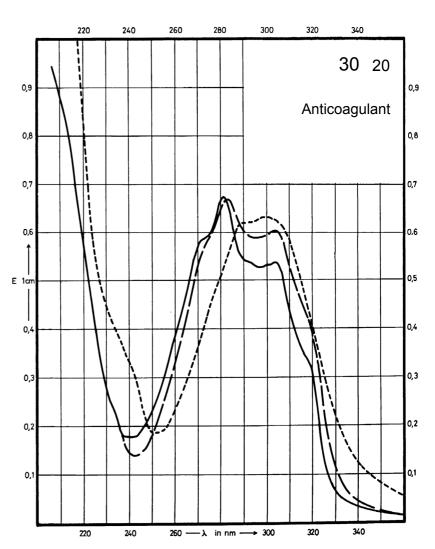

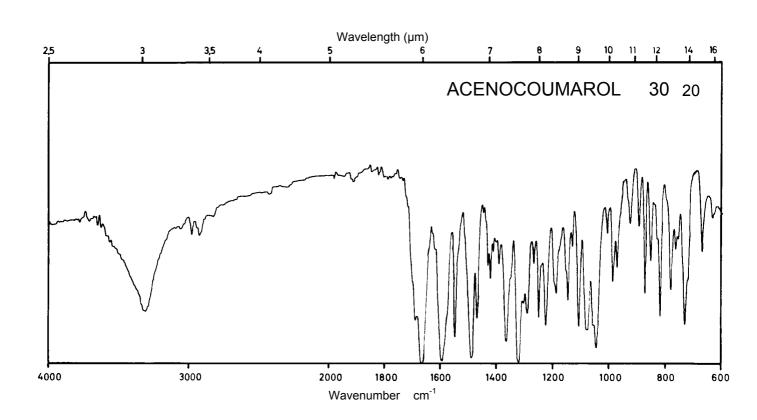





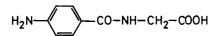

194.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm	272 nm	223 nm	273 nm
E ^{1%} 1cm	915	784	565	834
ε	17800	15200	11000	16200

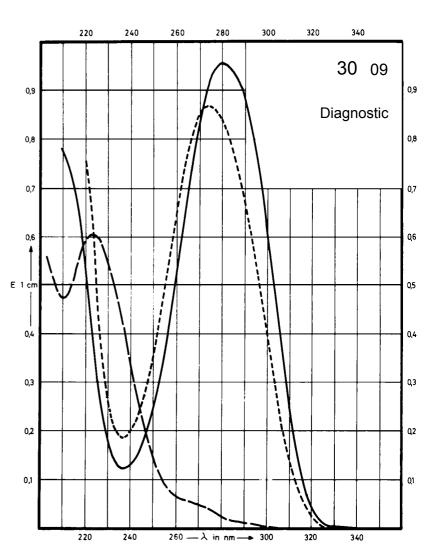

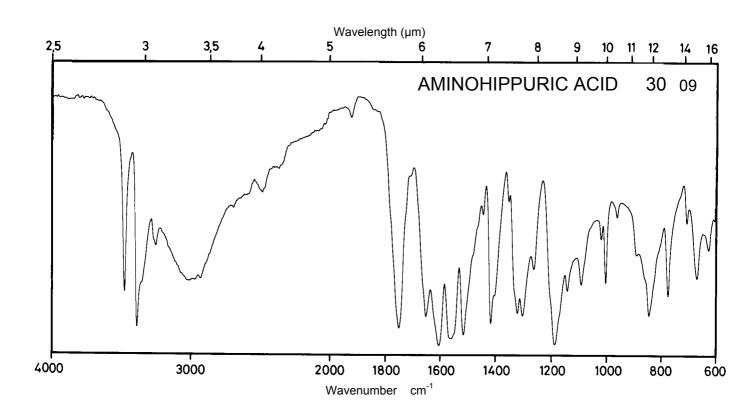



Name APOMORPHINE HYDROCHLORIDE

 M_{r}


303.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 275 nm		306 nm 272 nm	340 nm 264 nm
E ^{1%} 1cm	99 517		98 553	108 303
3	3010 15710		2980 16800	3280 9210

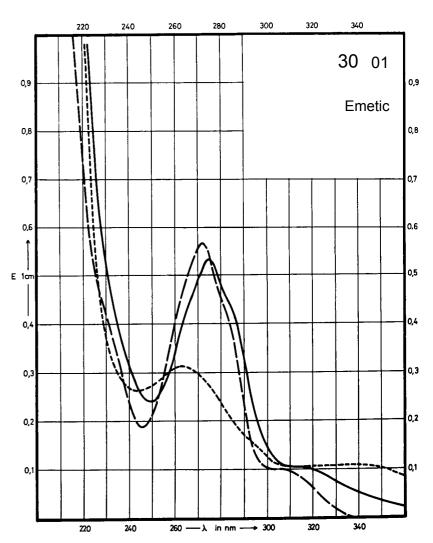

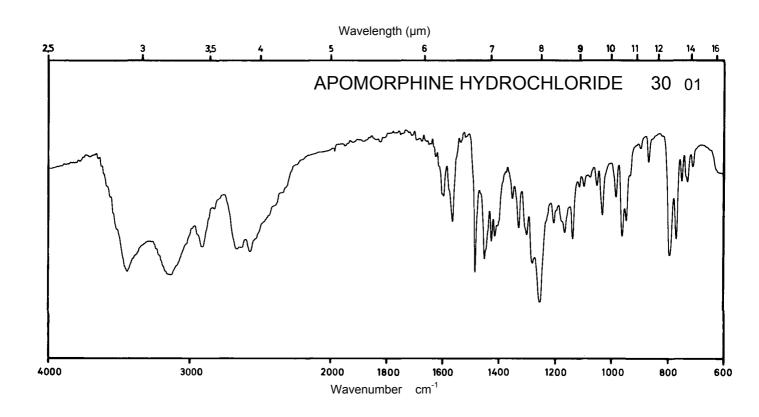

Name ANHYDROTETRA-CYCLINE HYDROCHLORIDE

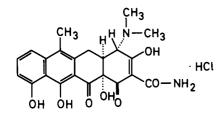


 M_{r}


462.9

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	430 nm 269 nm 224 nm		430 nm 274 nm 224 nm	430 nm 270 nm
E ^{1%} 1cm	223 1010 640		188 1110 629	249 900
ε	10330 46750 29670		8720 51400 29100	11530 41700

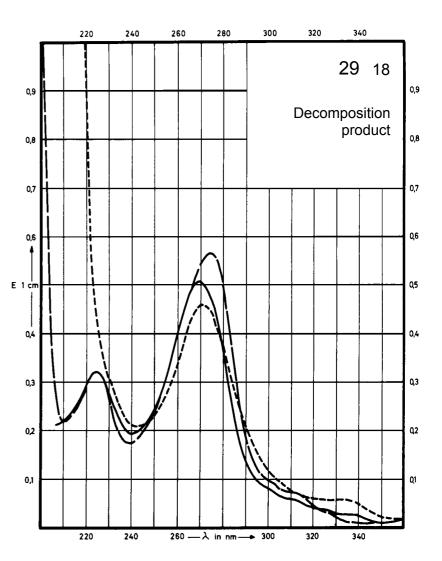

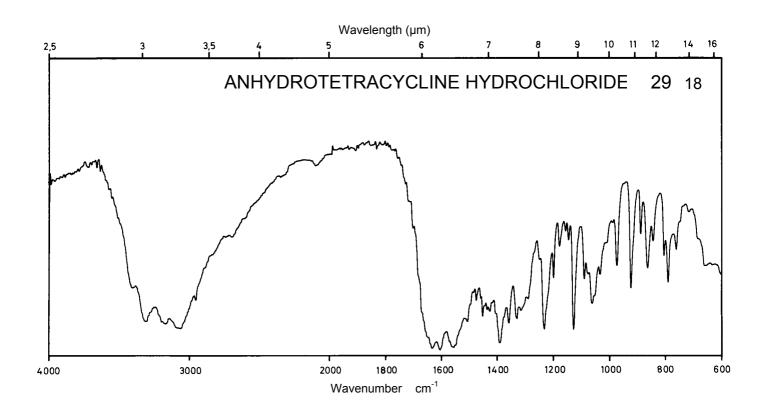





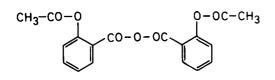

342.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 240 nm			
E ^{1%} 1cm	95 590			
ε	3250 20200			

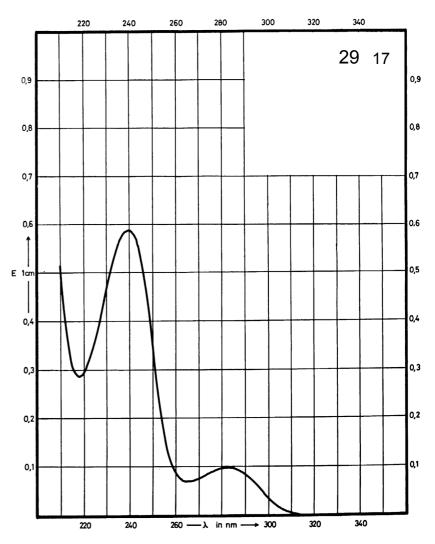

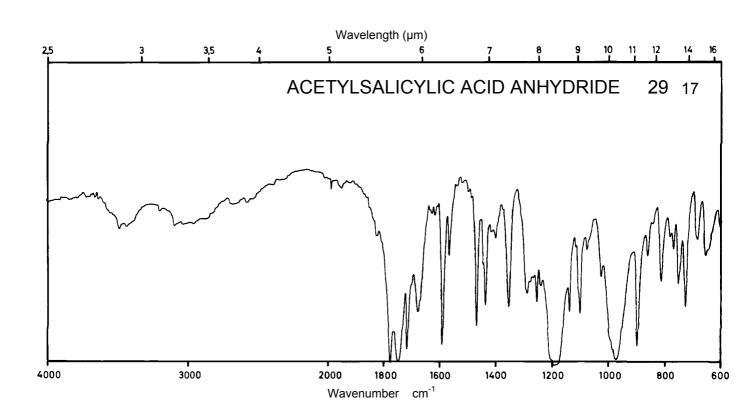


Name ACETAMINO PHENAZONE


$$H_3C$$
 $NH-CO-CH_3$
 H_3C-N
 N


 M_{r}


245.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		259 nm	260 nm
E ^{1%} 1cm	401		347	381
ε	9840		8510	9350

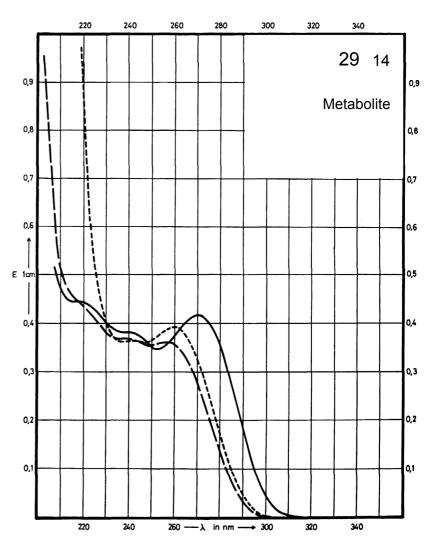

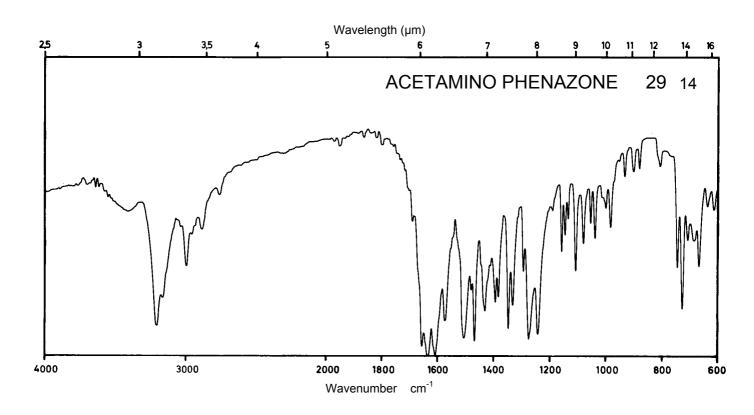

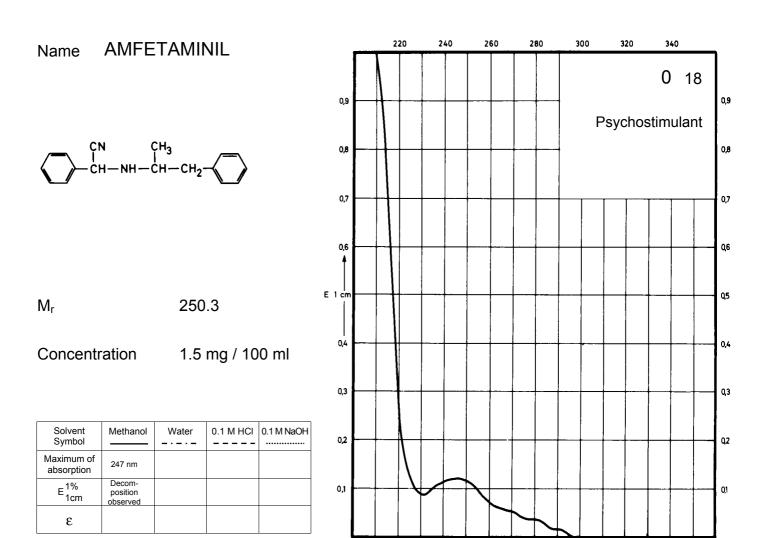

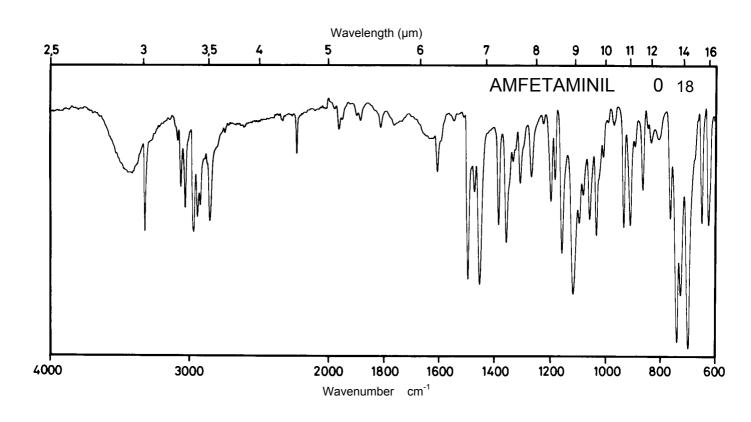

220

240

260 — λ in nm → 300

340


Name 4-AMINOPHENAZONE


H₃C NH₂

 M_{r}

203.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	251 nm		258 nm	244 nm
E ^{1%} 1cm	549		451	469
ε	11160		9170	9530

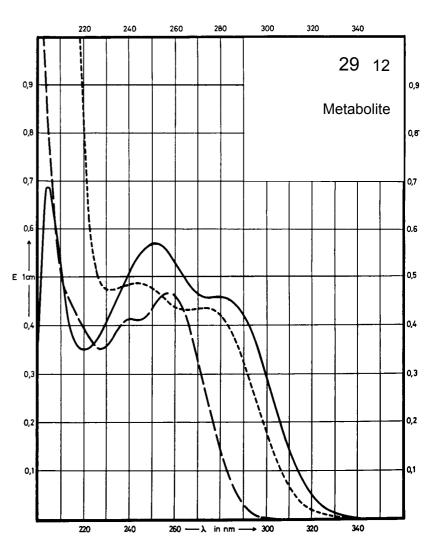

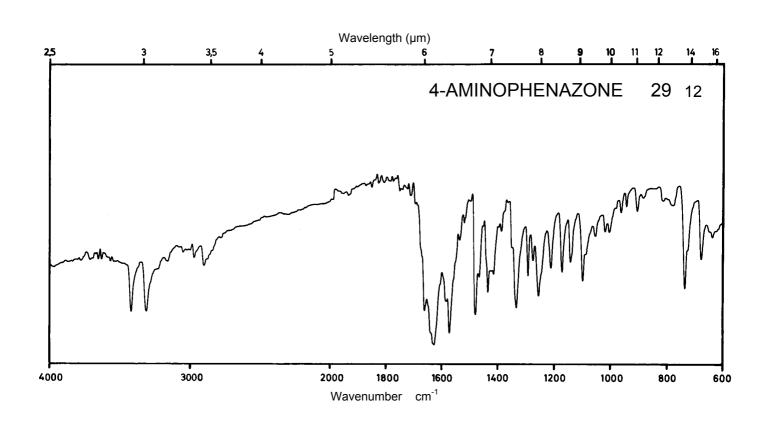

Name 4-AMINOPHENOL


OH OH


 M_{r}

109.1


Concentration 0.26 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	303 nm 233 nm		271 nm	266 nm
E ^{1%} 1cm	177 587		133	731
ε	1930 6400		1450	7980

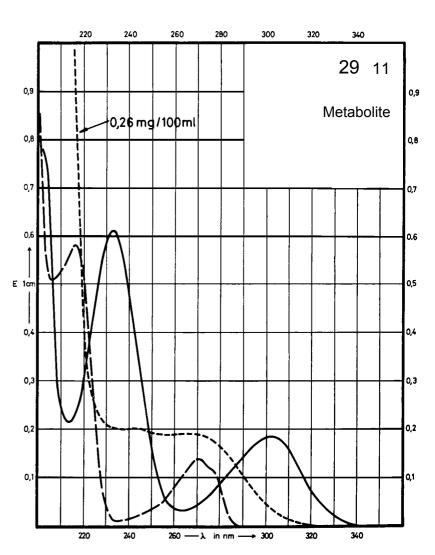

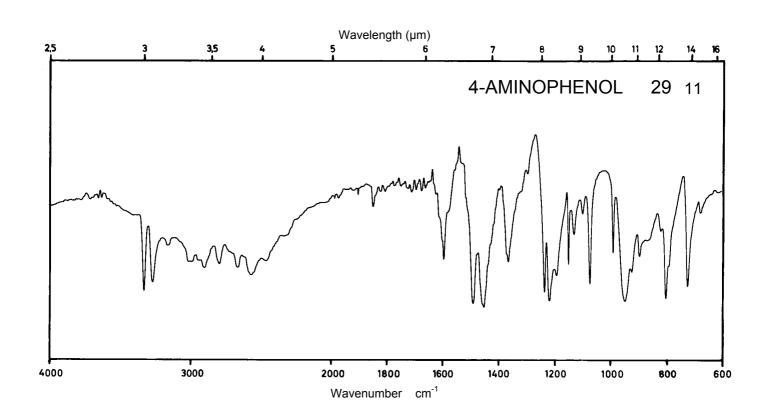




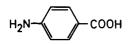


137.1


Concentration 0.5 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm		270 nm 225 nm	275 nm
E 1% 1cm	1335		95 854	1067
ε	18300		1300 11710	14630

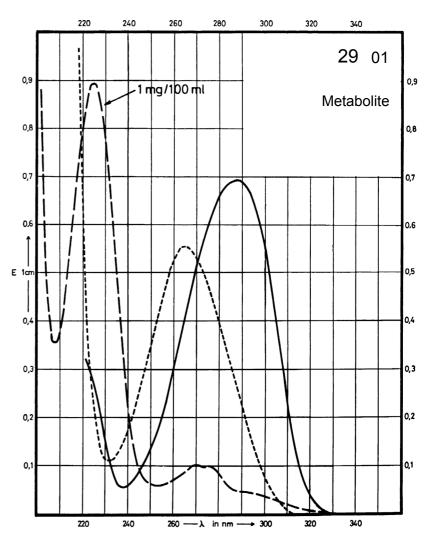

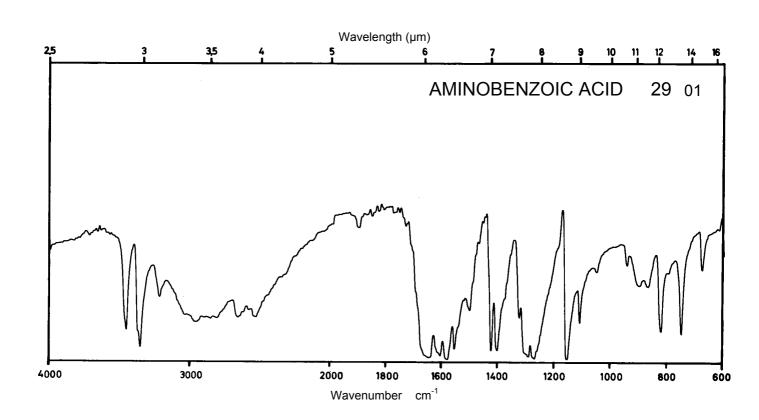



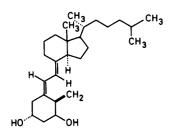



400.7

Concentration


0.95 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm			
E ^{1%} 1cm	430			
ε	17200			

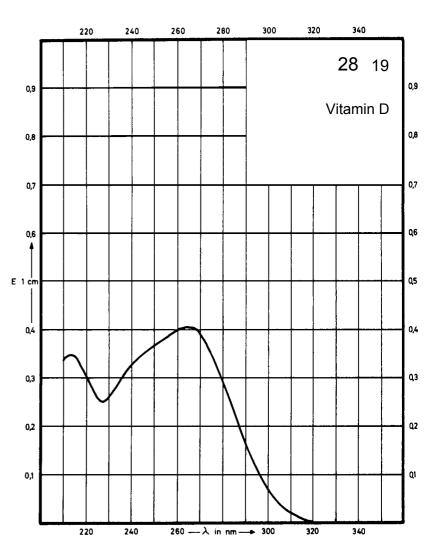

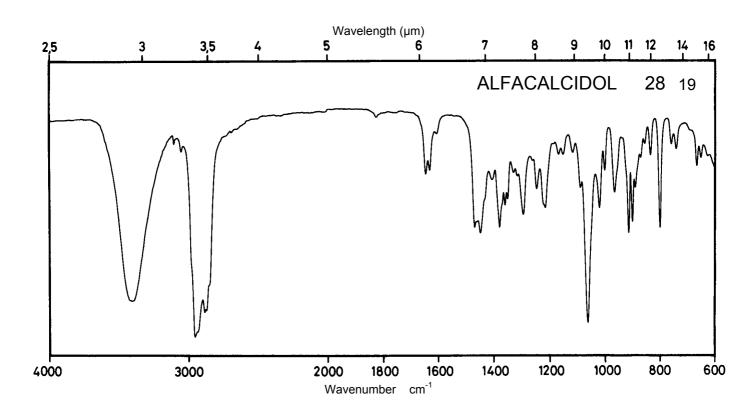





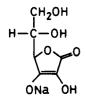


198.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			243 nm	
E ^{1%} 1cm			503	
ε			9960	

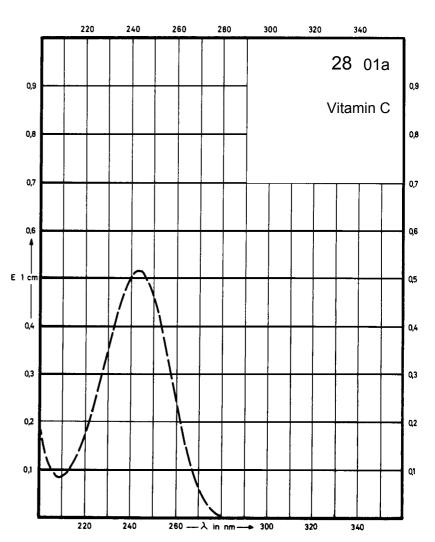

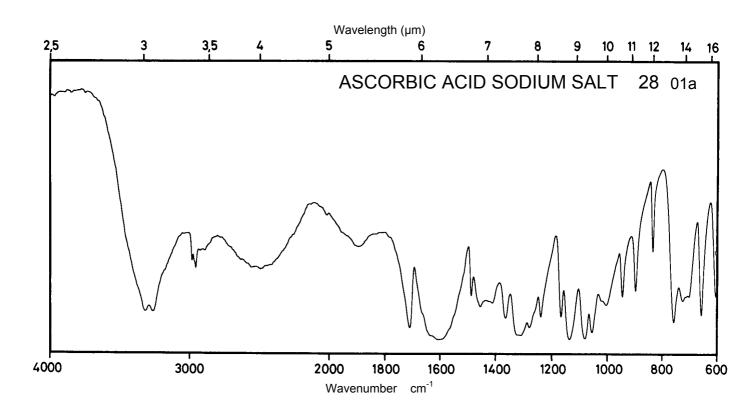


Name ASCORBIC ACID

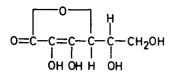


 M_{r}

176.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm		243 nm	
E ^{1%} 1cm	469		567	
3	8260		9980	

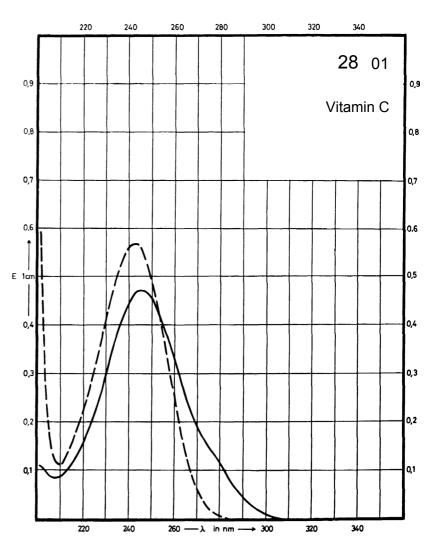

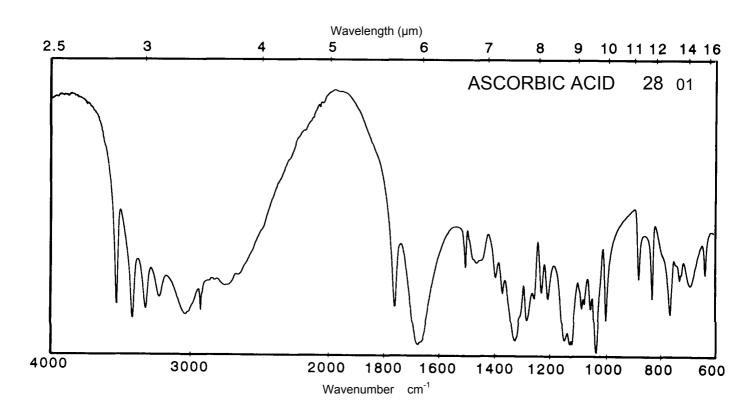

Name ANGIOTENSINAMIDE


 $\begin{array}{l} \operatorname{\mathsf{Asn}} - \operatorname{\mathsf{Arg}} - \operatorname{\mathsf{Val}} - \operatorname{\mathsf{Tyr}} - \operatorname{\mathsf{Val}} - \operatorname{\mathsf{His}} - \operatorname{\mathsf{Pro}} - \operatorname{\mathsf{Phe}} \\ \operatorname{\mathsf{NH}}_2 \end{array}$


 M_r

1031.2


Concentration 2.5 mg / 100 ml 50 mg / 100 ml

F				
Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		275 nm		
E ^{1%} 1cm		11.0		
3		1135		

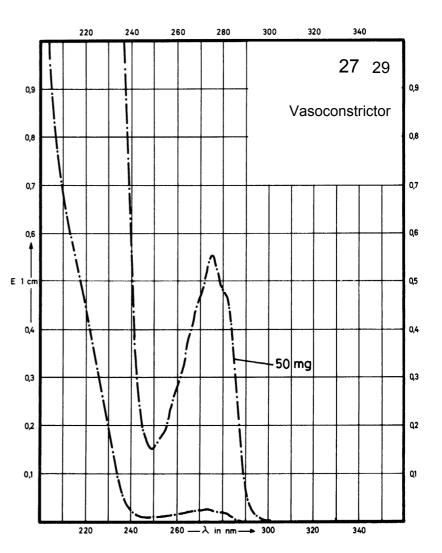

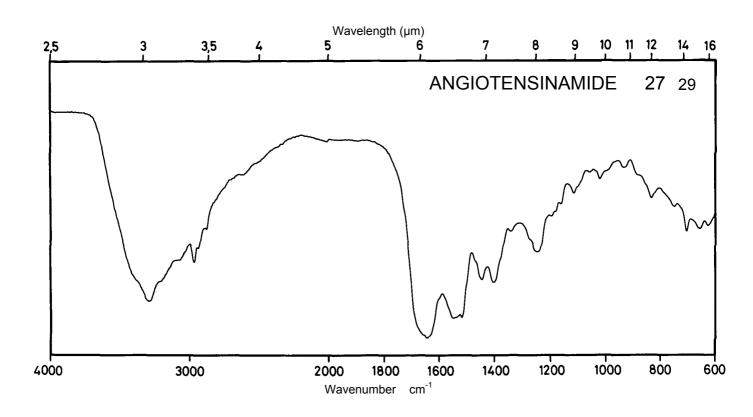





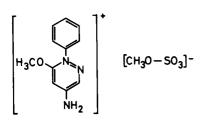


 M_r 313.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	290 nm 223 nm	288 nm 222 nm	288 nm 222 nm	Decom- position observed
E ^{1%} 1cm	479 430	468 462	468 462	
3	15000 13500	14700 14500	14700 14500	

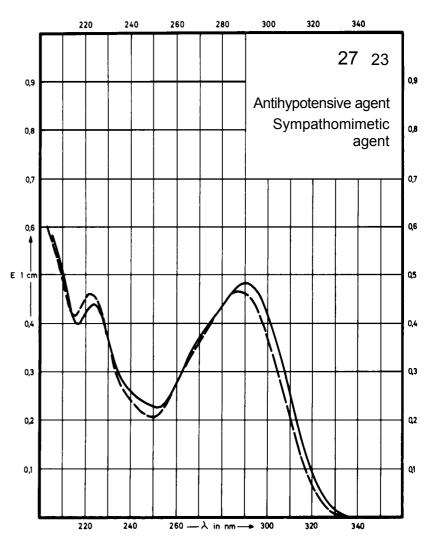

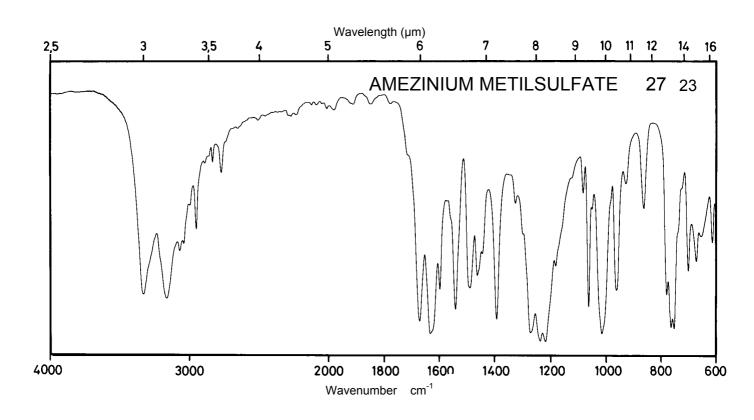





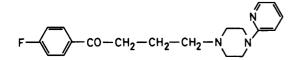


327.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	303 nm 247 nm		311 nm 242 nm	245 nm
E ^{1%} 1cm	123 863		264 744	723
ε	4000 28200		8650 24300	23700

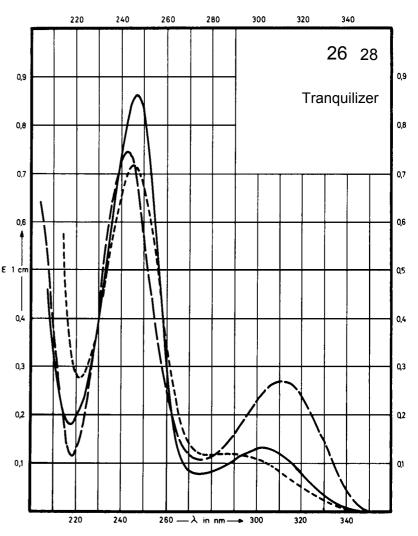

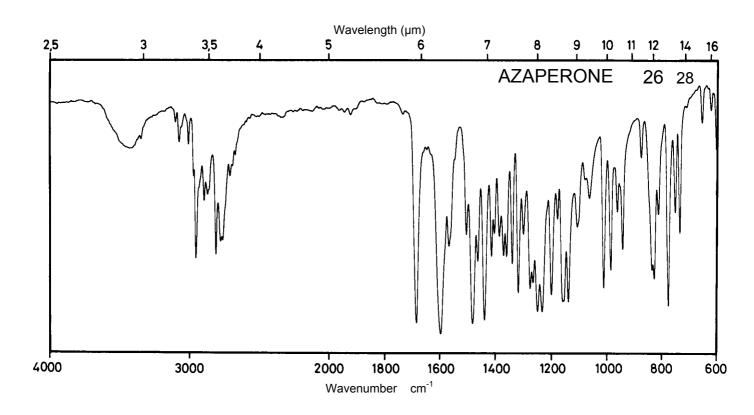



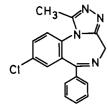



308.8

Concentration

0.56 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	221 nm		Decom- position observed	
E ^{1%} 1cm	1290			
ε	39800			

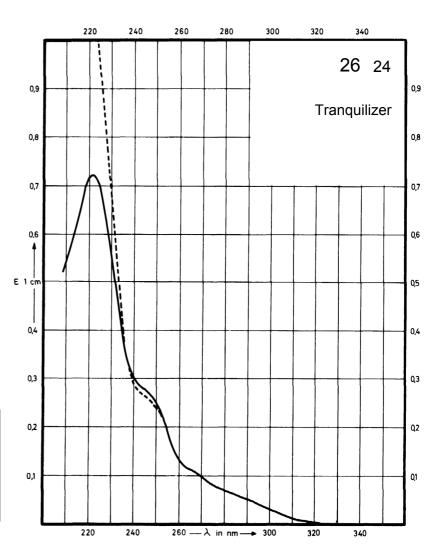

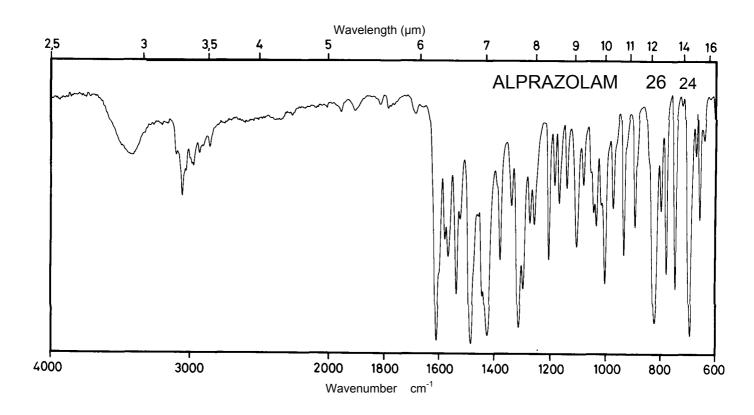





$$\begin{array}{c} \text{CH}_3 \\ \text{CH}_2 \\ \text{Br} - \overset{!}{\text{C}} - \text{CO} - \text{NH} - \text{CO} - \text{NH} - \text{CO} - \text{CH}_3 \\ \overset{!}{\text{CH}}_2 \\ \overset{!}{\text{CH}}_3 \end{array}$$

265.1

Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

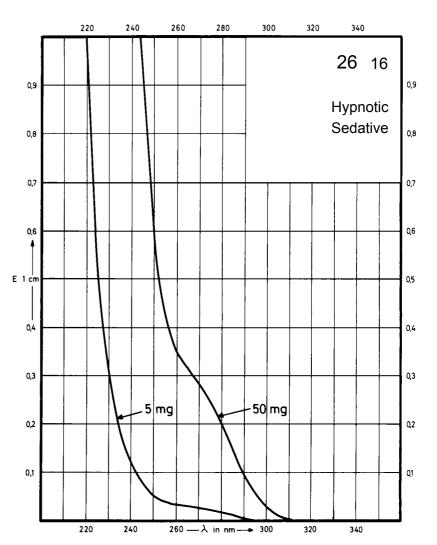

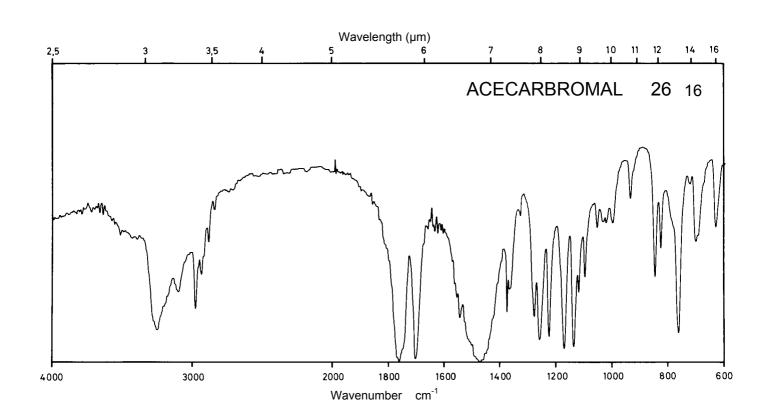







136.1

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	249 nm		250 nm	255 nm
E ^{1%} 1cm	566		555	521
ε	7700		7550	7090

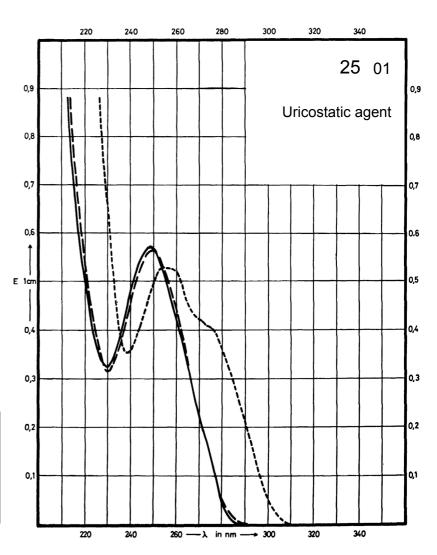

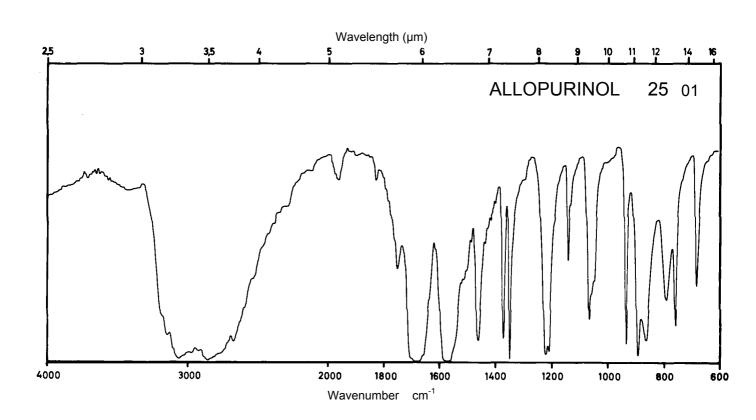





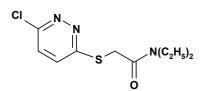

259.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm	255 nm	256 nm	256 nm
E ^{1%} 1cm	539	478	485	488
3	14000	12400	12600	12700

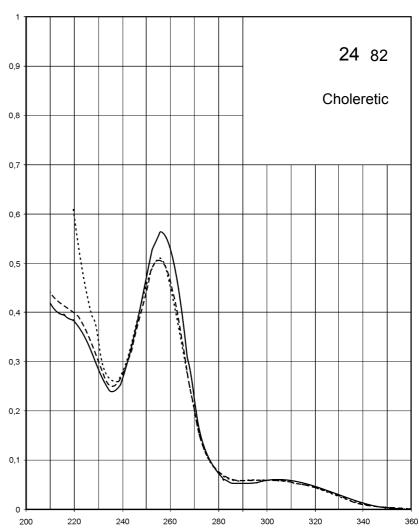

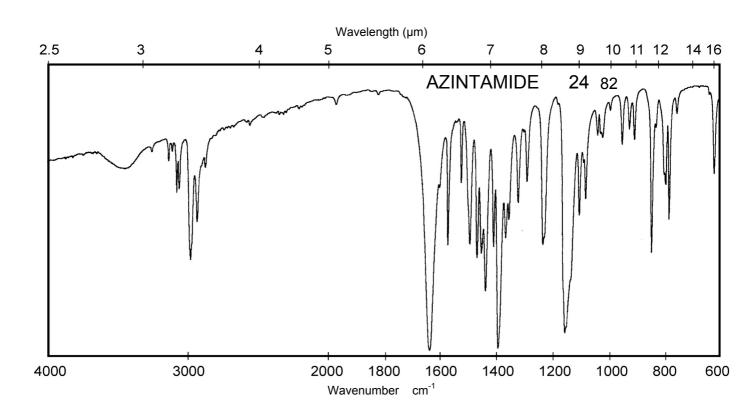



Name AMBUCETAMIDE HYDROCHLORIDE

 M_r


328.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm 237 nm		270 nm 236 nm	
E ^{1%} 1cm	39 389		35 380	
ε	1280 12800		1140 12500	

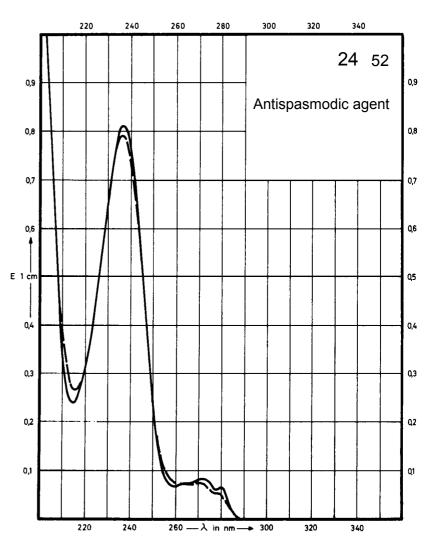

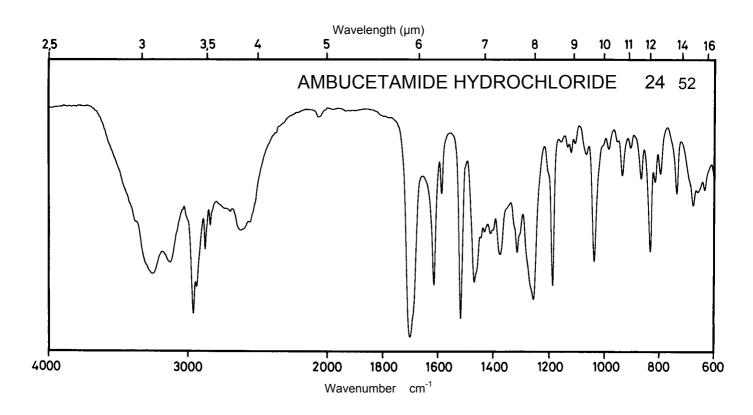





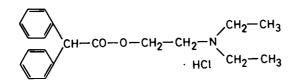


347.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 259 nm 253 nm	265 nm 259 nm 253 nm		
E ^{1%} 1cm	4.38 5.62 4.55	4.54 5.98 4.93		
ε	152 196	158 208 172		

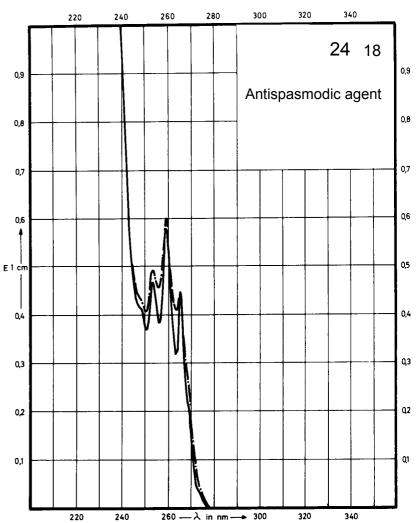

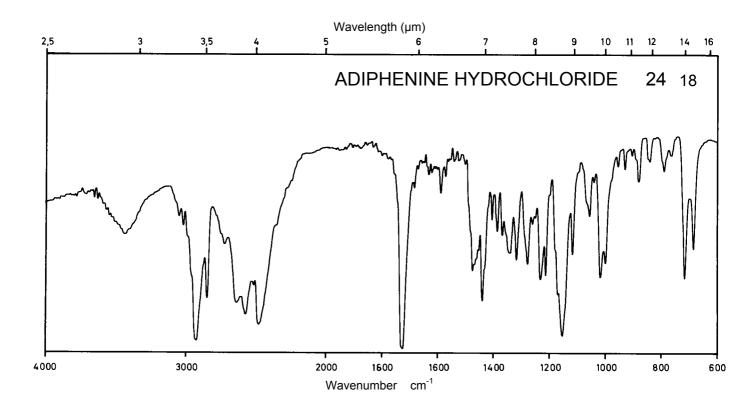





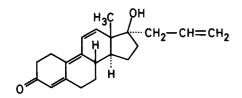


310.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	342 nm			
E ^{1%} 1cm	965			
ε	30000			

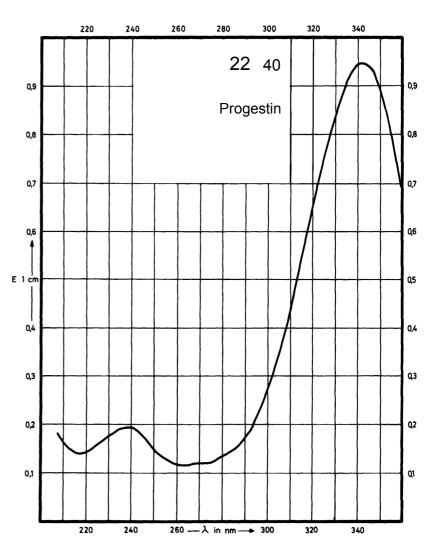

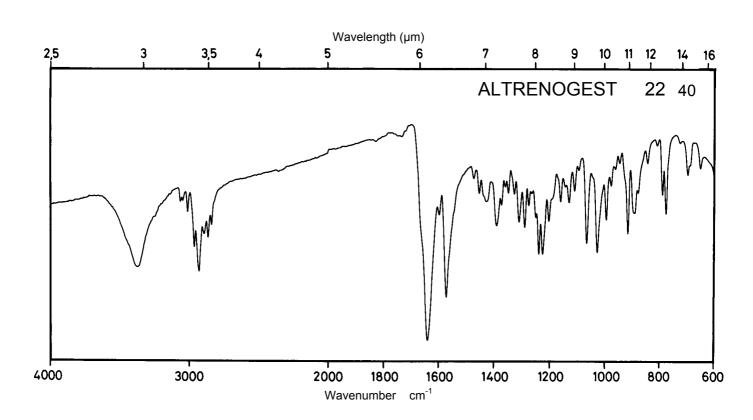





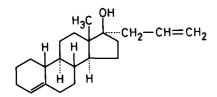


 M_r 300.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

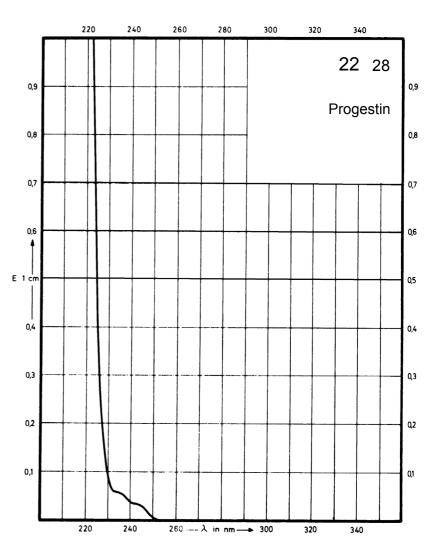

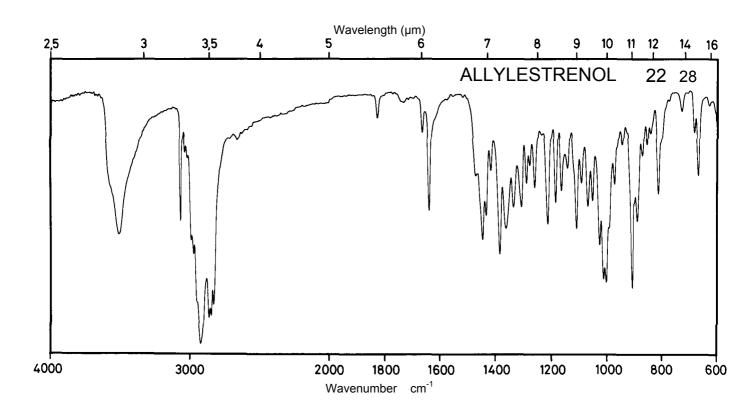





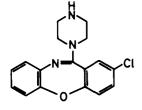

313.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	321 nm 298 nm 250 nm	298 nm	294 nm 252 nm	298 nm 250 nm
E ^{1%} 1cm	220 224 396	256	311 355	266 375
ε	6900 7030 12430	8020	9750 11150	8350 11770

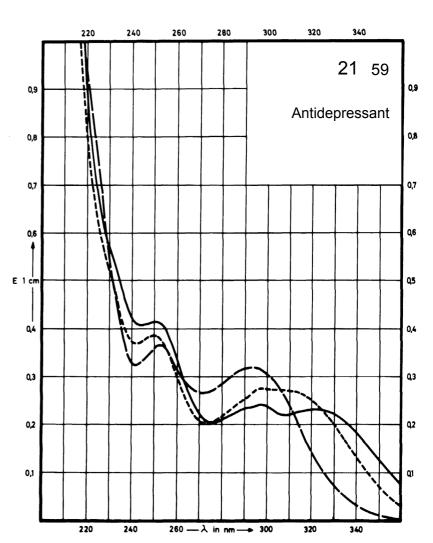


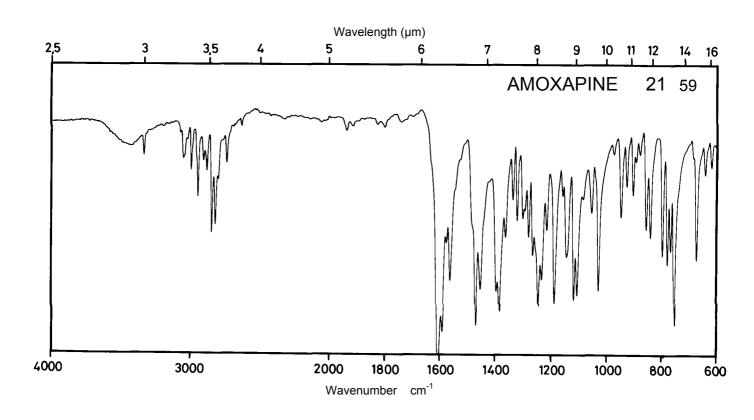


Name AMITRIPTYLINOXIDE

 M_r 329.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm	238 nm	238 nm	238 nm
E ^{1%} 1cm	441	421	420	426
ε	14500	13900	13800	14000

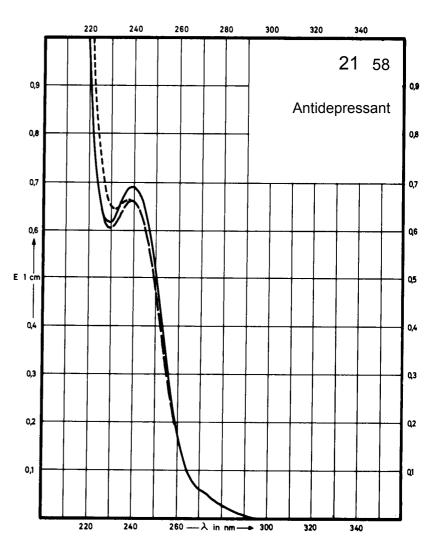

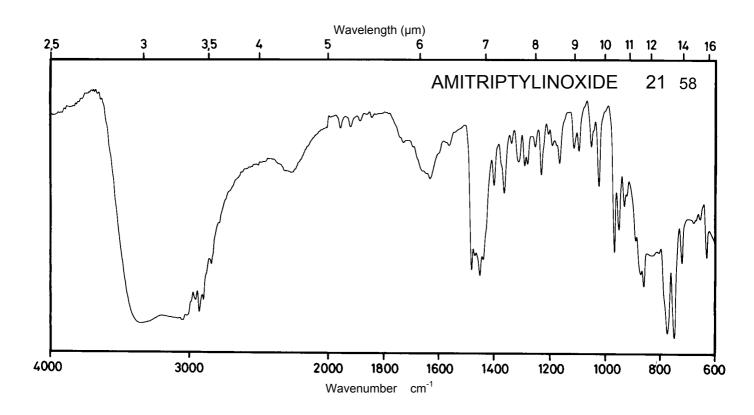


Name AMFETAMINE SULFATE


$$\left[\begin{array}{c} -\text{CH}_2 - \text{CH} - \text{CH}_3 \\ + \text{NH}_3 \end{array}\right]_2 \cdot \text{SO}_4^{--}$$

 M_{r}


368.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm	263 nm 257 nm 252 nm	263 nm 257 nm 252 nm	
E ^{1%} 1cm	7.5 9.7 7.6	7.5 9.9 7.9	7.5 9.9 7.9	
ε	278 359 281	278 365 291	278 365 291	

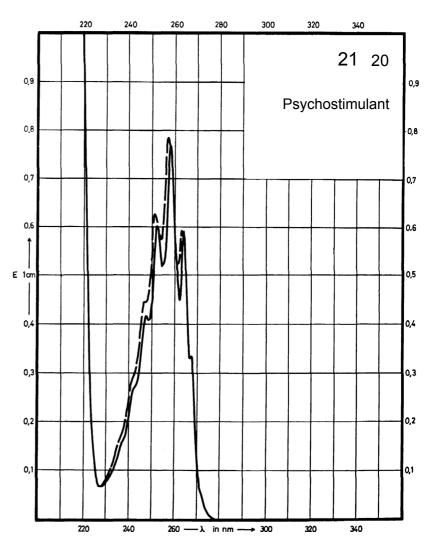

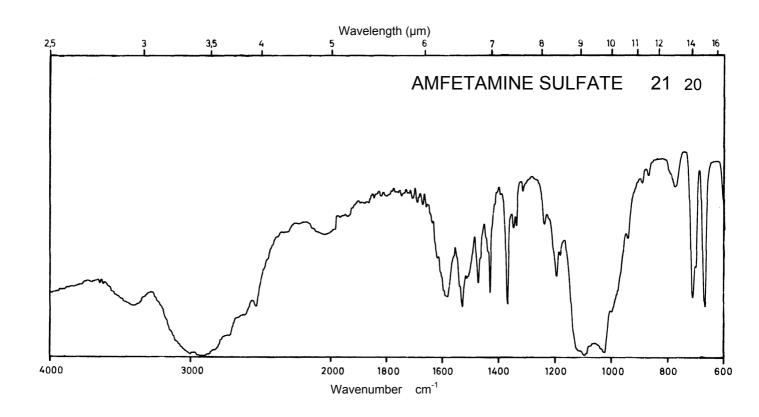


Name AMITRIPTYLINE HYDROCHLORIDE

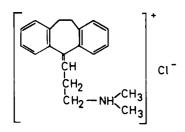



 M_{r}


313.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm		239 nm	237 nm
E 1% 1cm	464		441	388
ε	14560		13840	12180

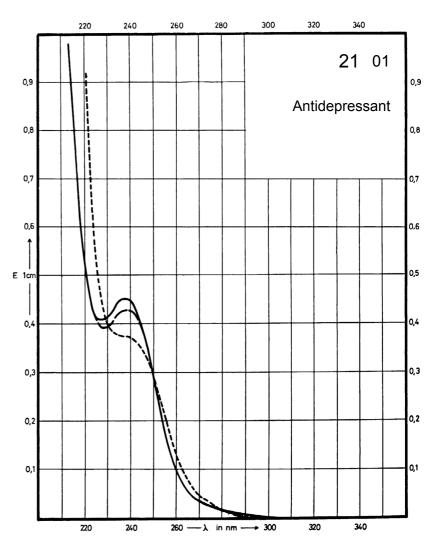


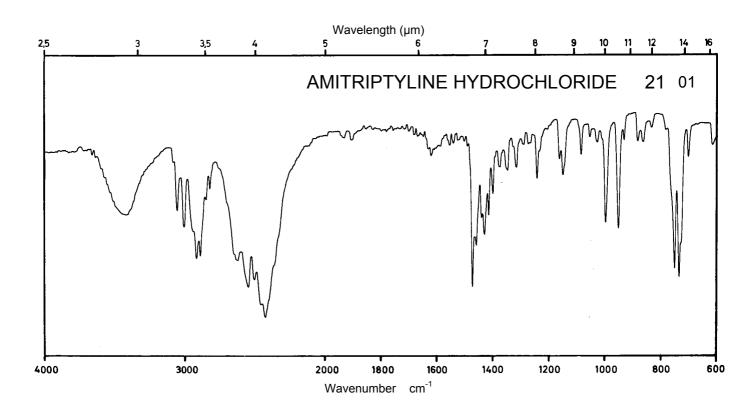

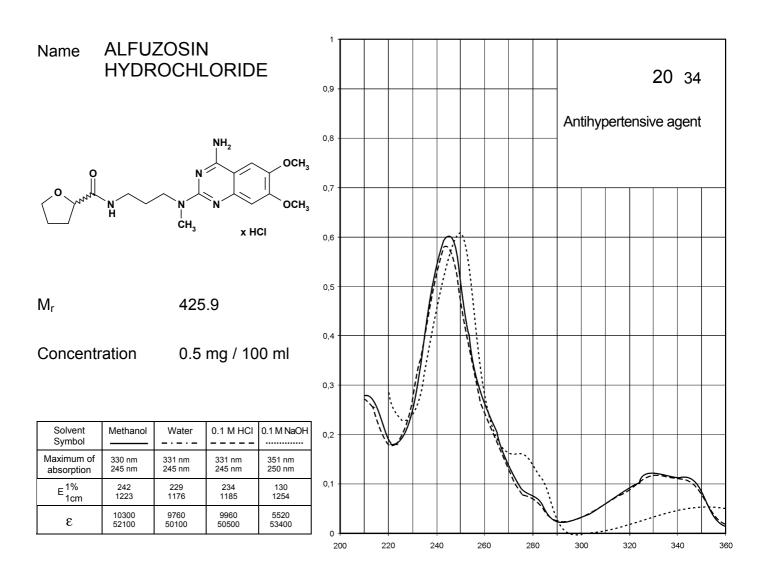

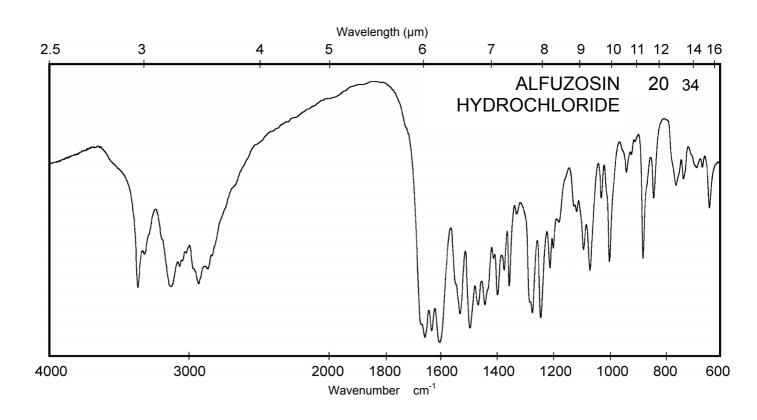



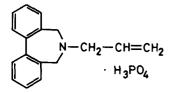



333.3

Concentration 1.5 r


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	247 nm		247 nm	242 nm
E ^{1%} 1cm	436		463	379
ε	14500		15400	12600

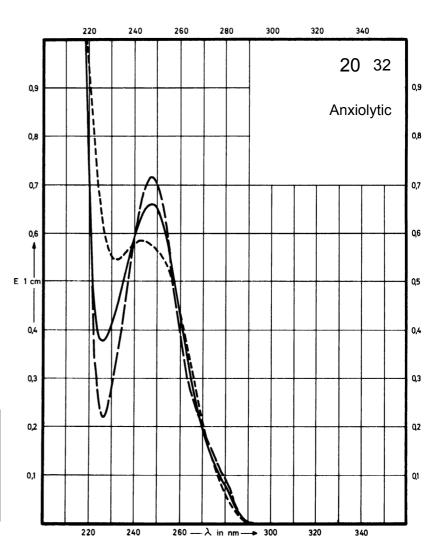

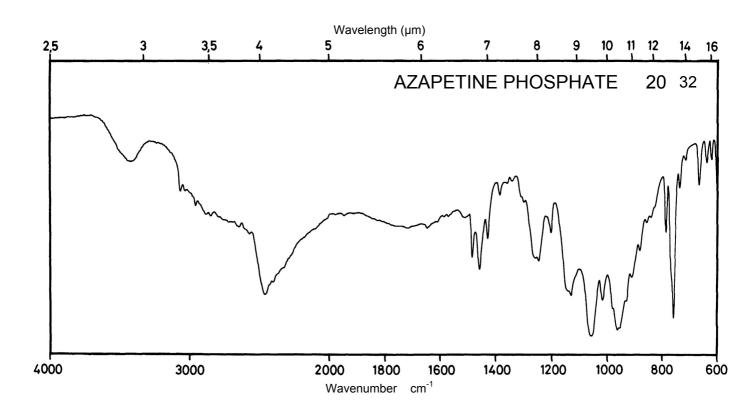




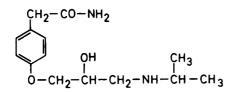


 M_r 266.3

Concentration 2.4 mg / 100 ml 12 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 227 nm		274 nm 225 nm	275 nm
E ^{1%} 1cm	54 375		47.5 349	49
3	1430 10000		1260 9290	1300

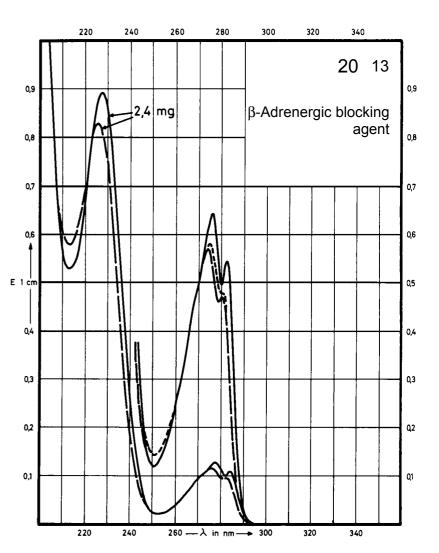

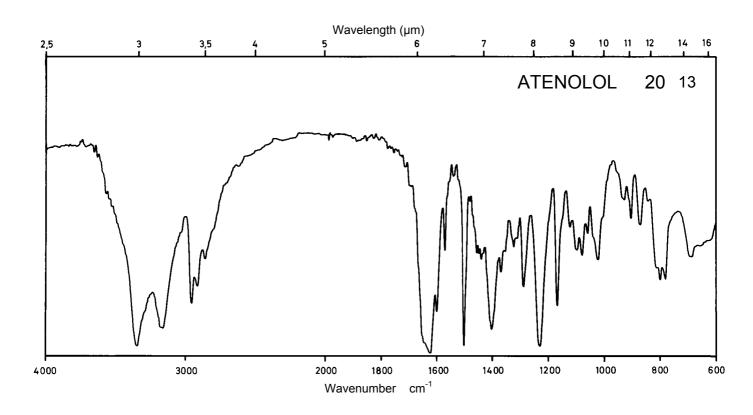





$$\begin{array}{c} \text{CH}_{3}\text{-CH} - \text{CH} \\ \text{NH} \\ \text{CO-CH}_{3} \quad \text{CH}_{2} \\ \text{CH}_{3}\text{-CH}_{2}\text{-CH}_{2}\text{-CO-NH} - \begin{array}{c} \text{CO-CH}_{2}\text{-CH} \\ \text{-} \text{HCI} \quad \text{OH} \end{array}$$

 M_{r} 372.9


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 328 nm 235 nm 320 nm 234 nm 320 nm 234 nm absorption E^{1%} 1cm 69 708 68 583 68 591 2580 26410 2540 22040 ε

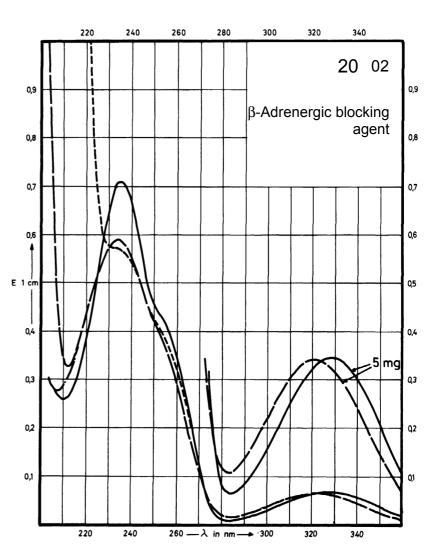

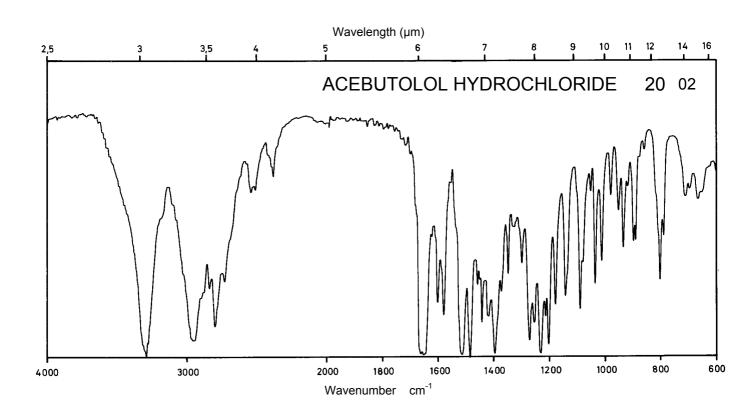




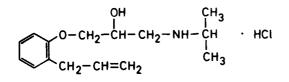


 M_r


285.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 271 nm		276 nm 270 nm	271 nm
E ^{1%} 1cm	58 64		53 60	62
ε	1650 1840		1520 1710	1770

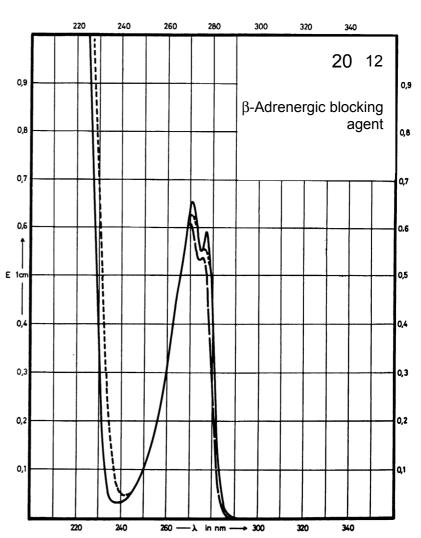

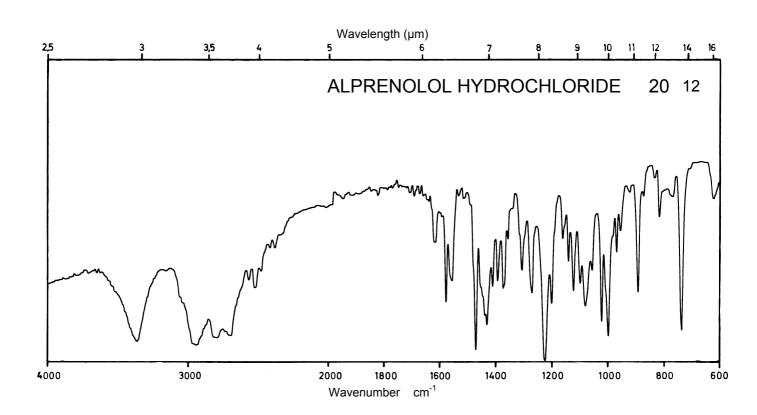





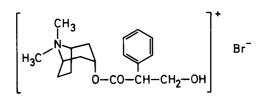


384.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm		257 nm	257 nm
E ^{1%} 1cm	5.0		5.2	5.4
3	190		198	209

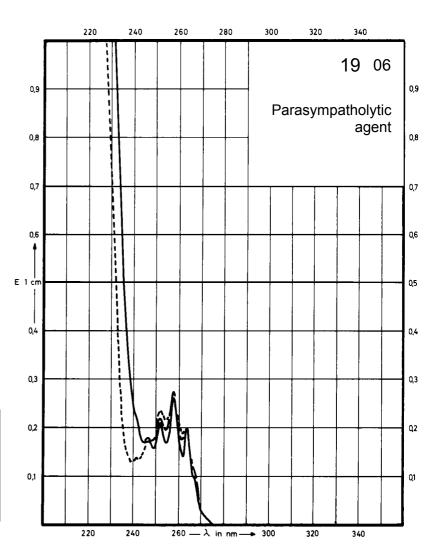

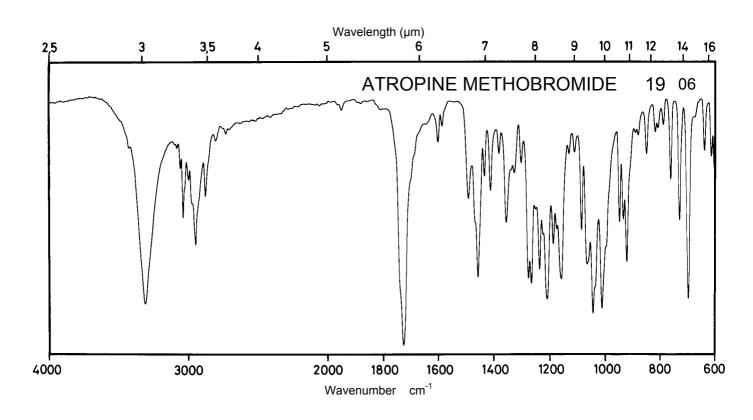





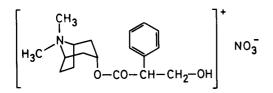

 M_{r} 366.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm			
E ^{1%} 1cm	5.2			
ε	190			

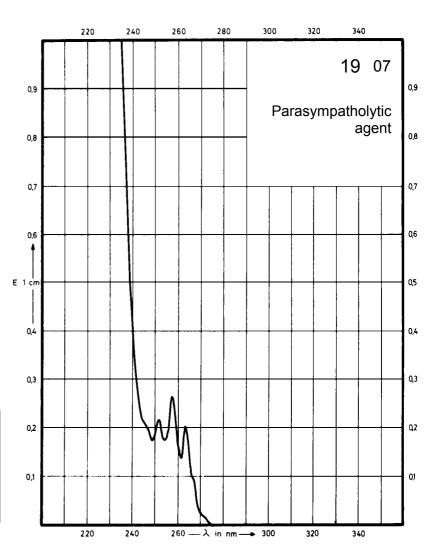

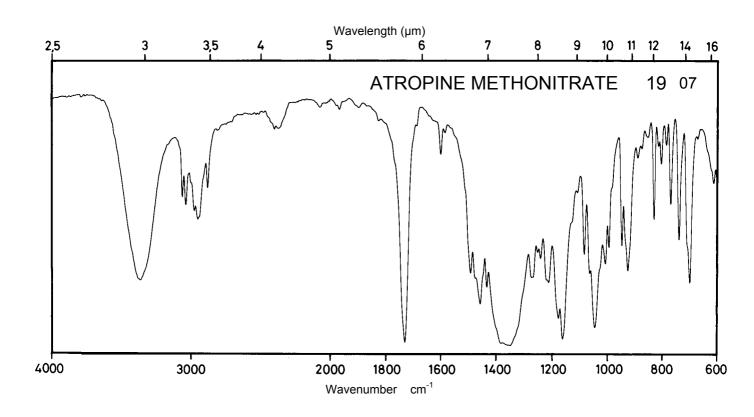



Name ATROPINE SULFATE


 M_{r}


694.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm		257 nm	257 nm
E ^{1%} 1cm	5.6		5.7	7.3
ε	390		400	510

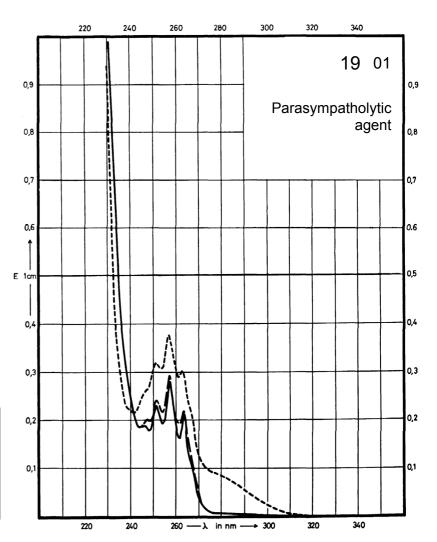

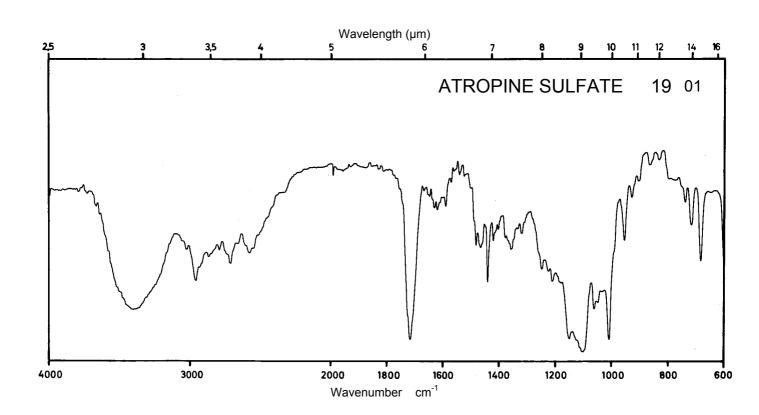

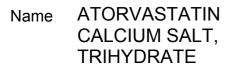


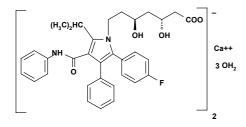


 M_{r}

1209.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	247 nm	241 nm	242 nm	240 nm
E ^{1%} 1cm	429	377	369	377
ε	51800	45600	44600	45600

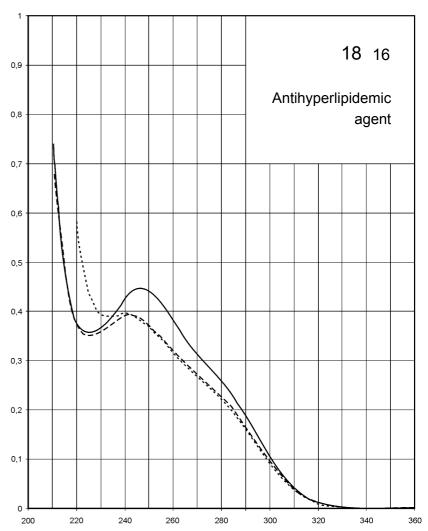

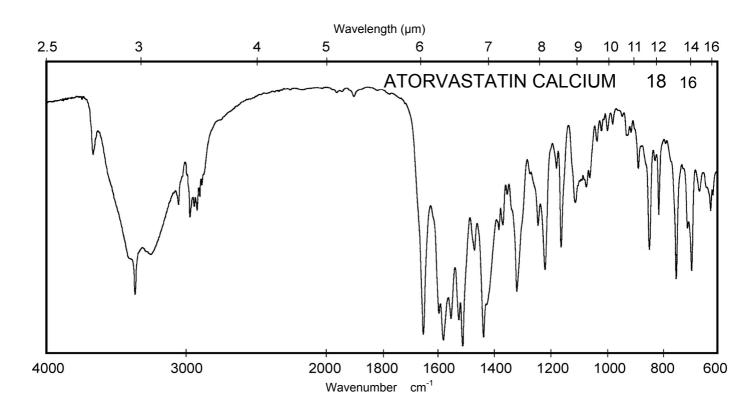



Name ACIPIMOX

H₃C N COOH


 M_r 154.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	232 nm 273 nm		230 nm 268 nm	224 nm 264 nm
E ^{1%} 1cm	1330 610		1452 598	1330 637
ε	20500 9400		22370 9220	20500 9810

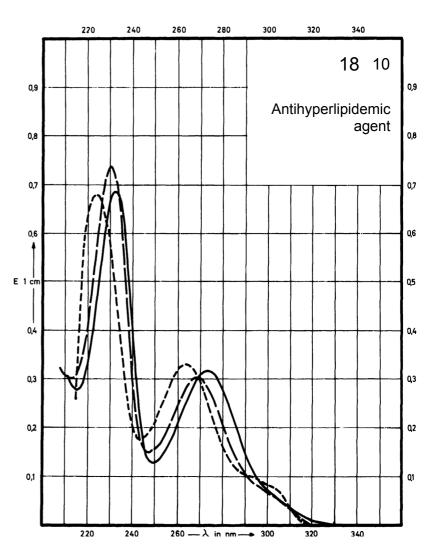

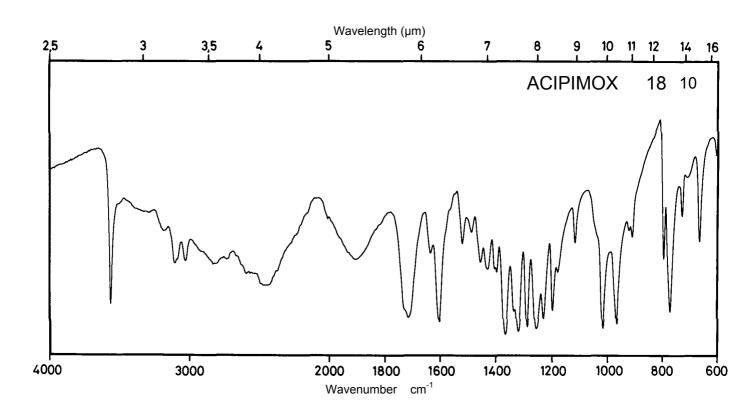





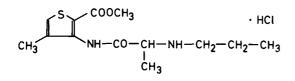

 M_r 320.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm		270 nm	
E ^{1%} _{1cm}	372		305	
ε	11940		9790	

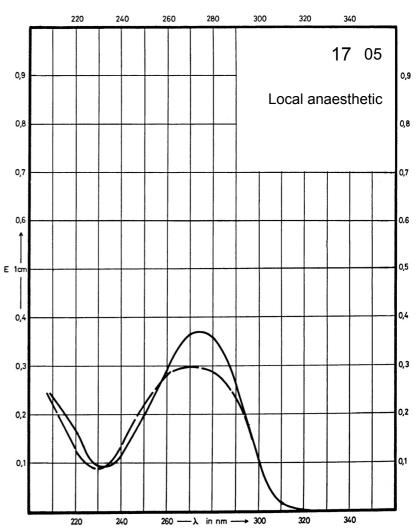

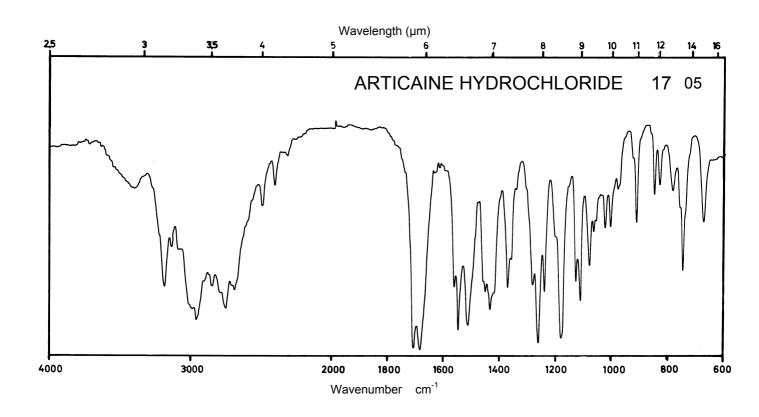



Name ALLOBARBITAL

$$OC$$
 NH
 CO
 CH_2
 CH_2

 M_r 208.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				250 nm
E ^{1%} 1cm				303
3				6300

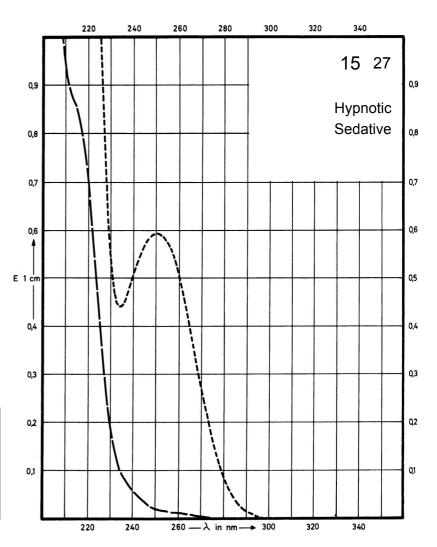

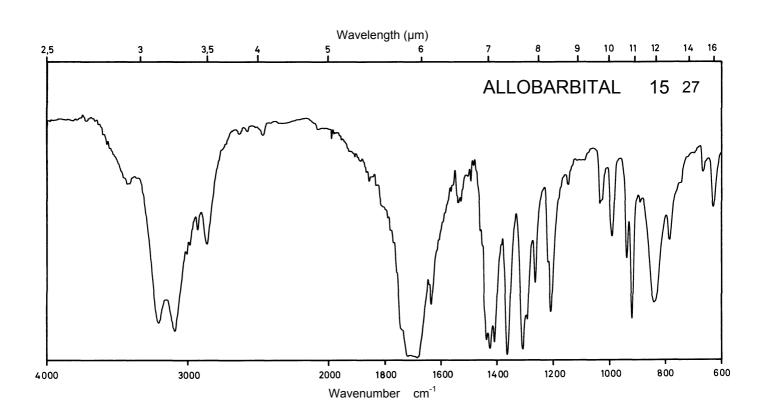


Name APROBARBITAL


CH₂=CH-CH₂ CO-NH CH₃-CH CO-NH CH₃

 M_{r}


210.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				243 nm
E ^{1%} 1cm				333
ε				7010

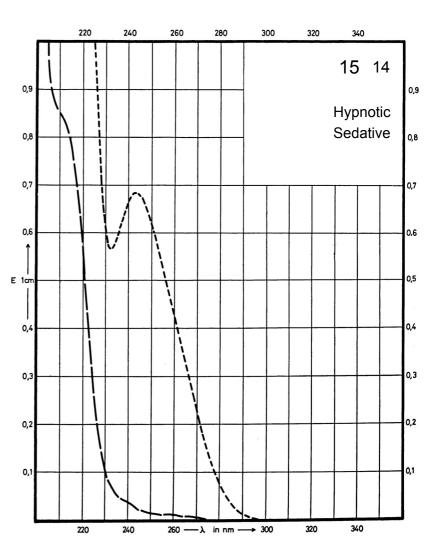

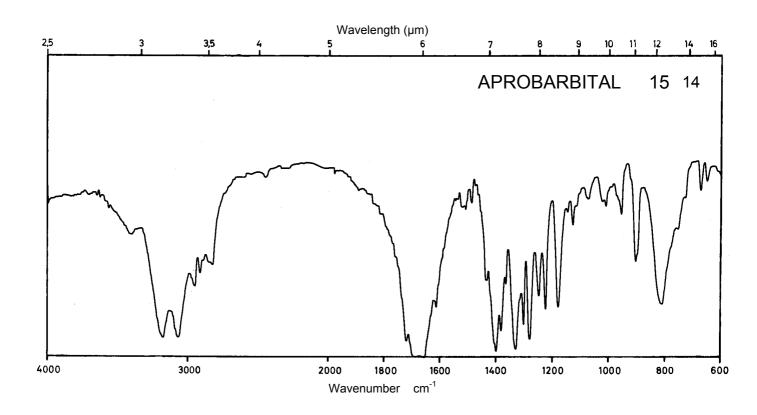



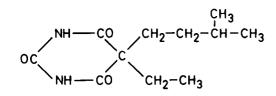



226.3

Concentration

2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				244 nm
E ^{1%} 1cm				318
ε				7200

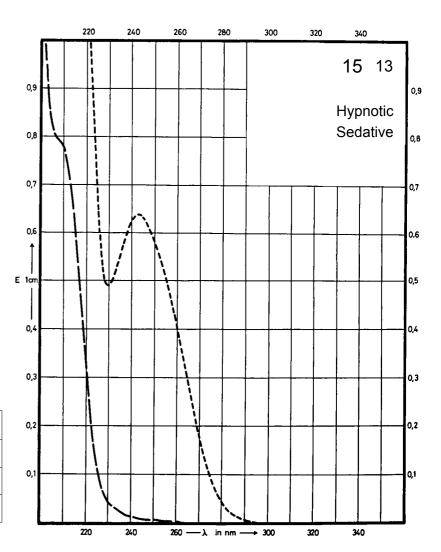

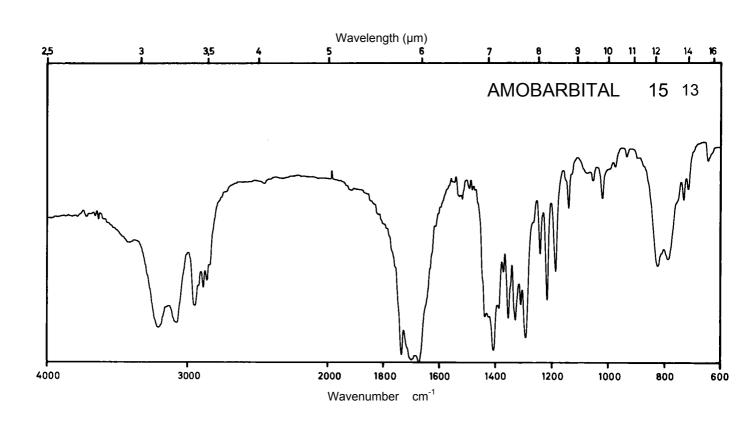



Name AMLODIPINE BESILATE

 M_{r}


567.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	361 nm	360 nm	360 nm	360 nm
	238 nm	239 nm	239 nm	239 nm
E ^{1%}	119	117	118	112
1cm	342	319	319	309
3	6770	6650	6670	6360
	19370	18080	18080	17500

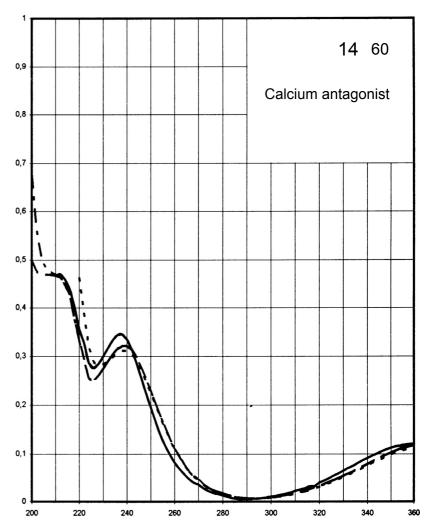

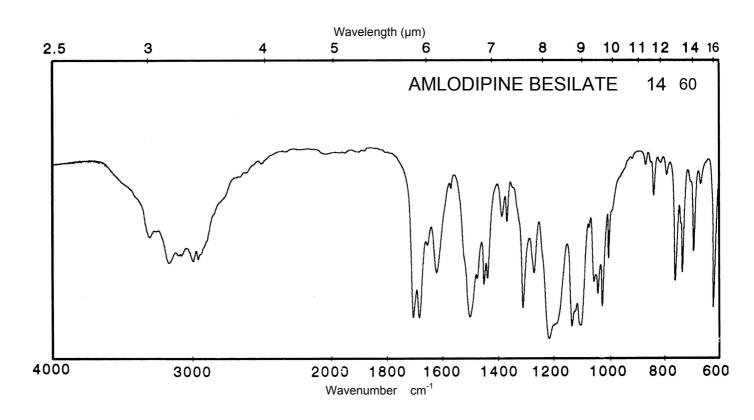





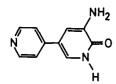


187.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	326 nm 244 nm		317 nm 273 nm	330 nm 234 nm
E ^{1%} 1cm	807 1180		1080 445	913 955
ε	15100 22100		20200 8300	17100 17900

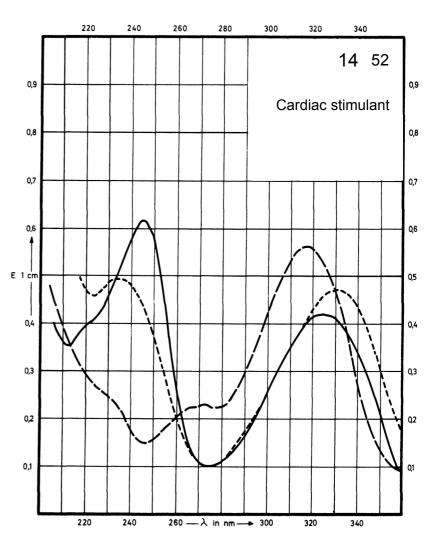

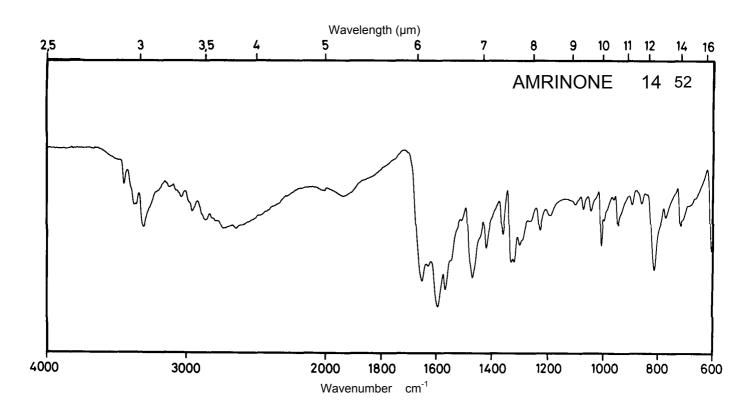




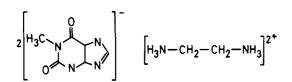


 M_r 420.4


Concentration 1.3 mg / 100 ml (anhydrous)

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm		270 nm	275 nm
E ^{1%} _{1cm}	535		543	648
ε	22500		22800	27200

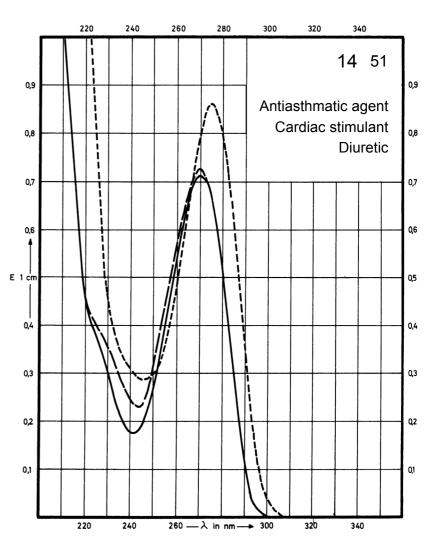

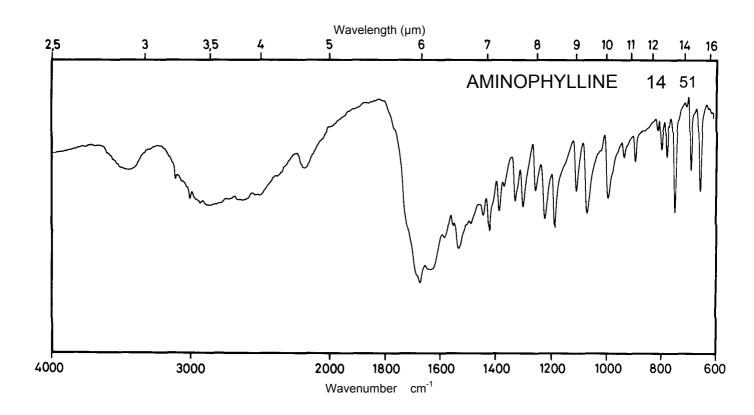




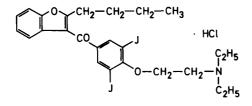


 M_r


681.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm		Decom- position observed	
E ^{1%} 1cm	553			
ε	37700			

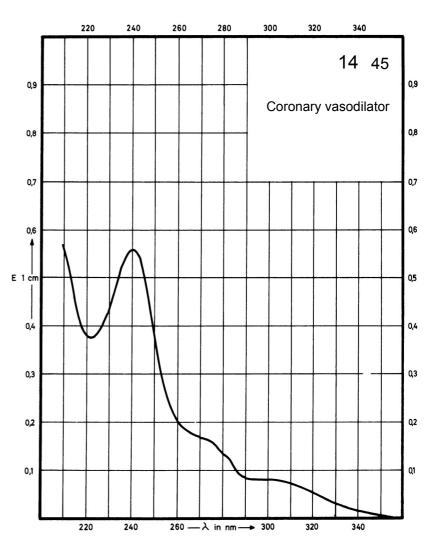

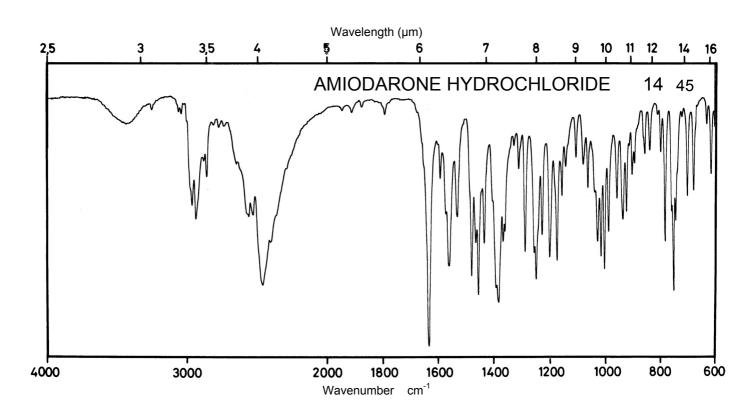


Name APRINDINE HYDROCHLORIDE

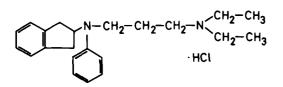


 M_{r}


359.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm 296 nm	253 nm	271 nm 264 nm	
E ^{1%} 1cm	349 40	211	32 34	
3	12530 1420	7560	1140 1210	

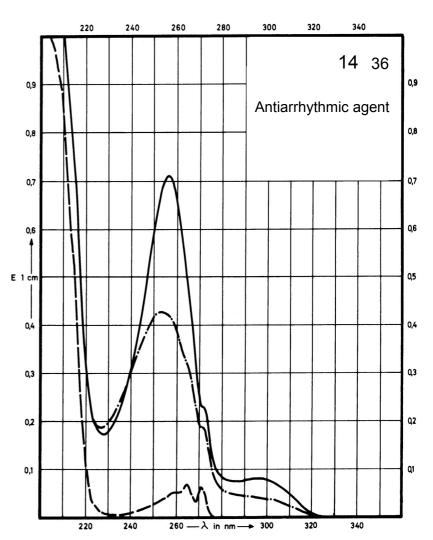

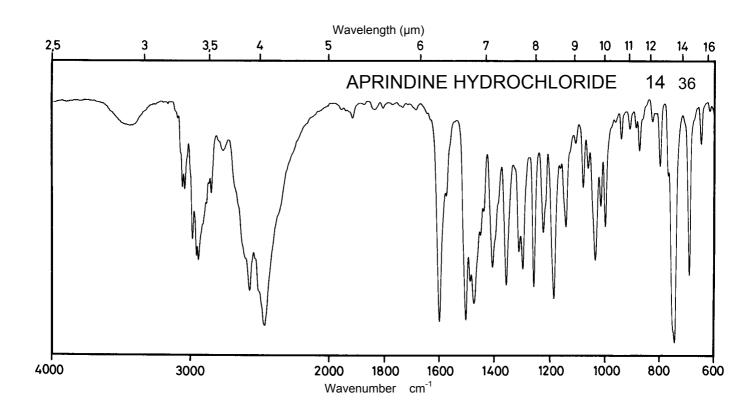



Name ACETYLDIGOXIN, β -ISOMER

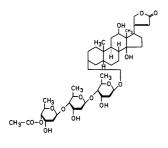


 M_r 823.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	219 nm	221 nm		
E ^{1%} 1cm	186	183		
3	15300	15050		

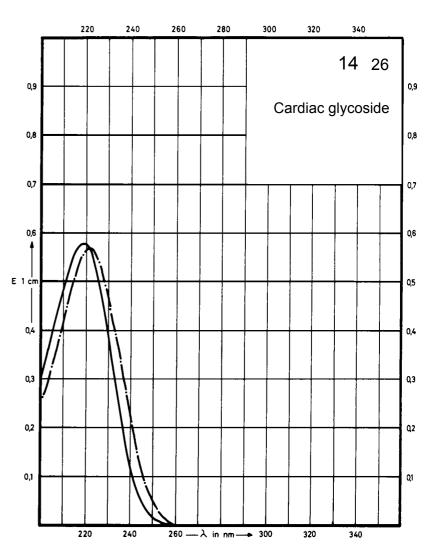

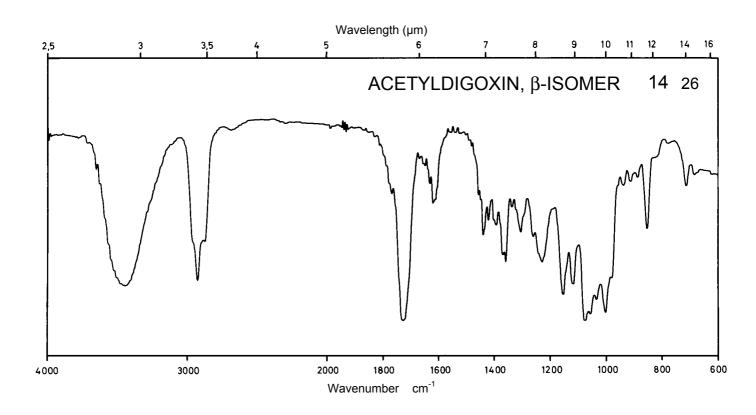





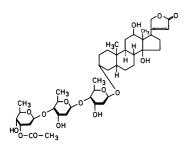


823.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	219 nm	221 nm		
E ^{1%} 1cm	186	183		
ε	15300	15050		

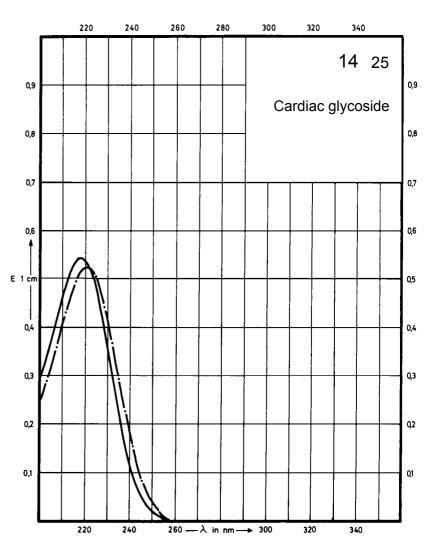

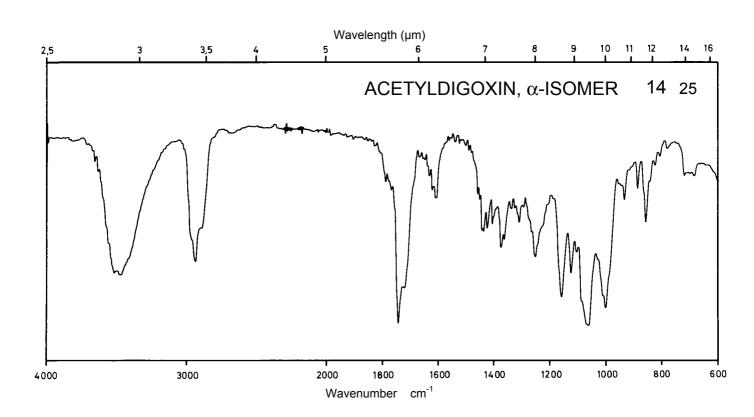

Name AJMALINE

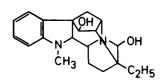


 M_{r}


326.4

Concentration


2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	292 nm 249 nm		289 nm 245 nm	290 nm 247 nm
E ^{1%} 1cm	89 262		75 221	76 220
ε	2890 8560		2440 7200	2480 7170

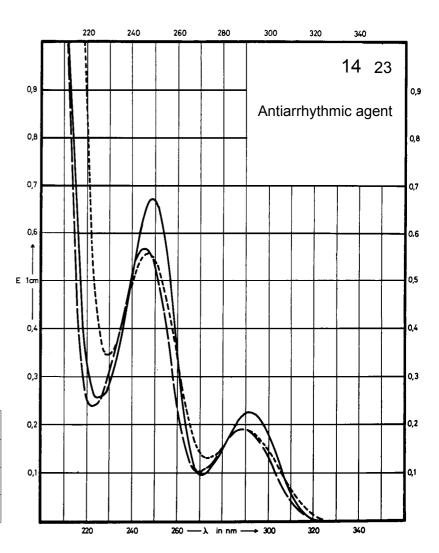

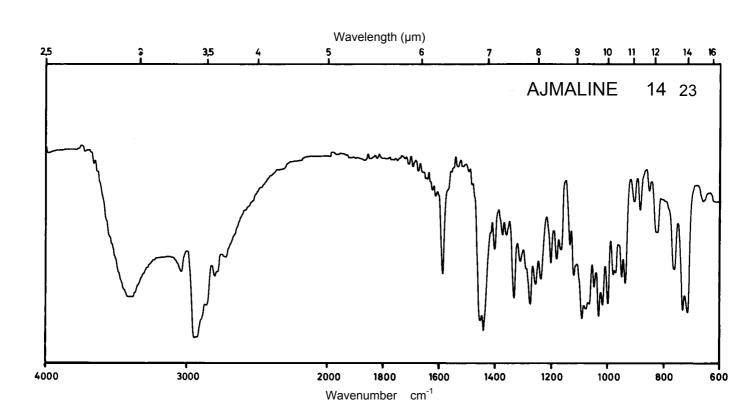


Name ADENOSINE

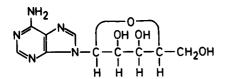


 M_{r}

267.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm		256 nm	260 nm
E ^{1%} 1cm	566		540	555
ε	15120		14430	14830

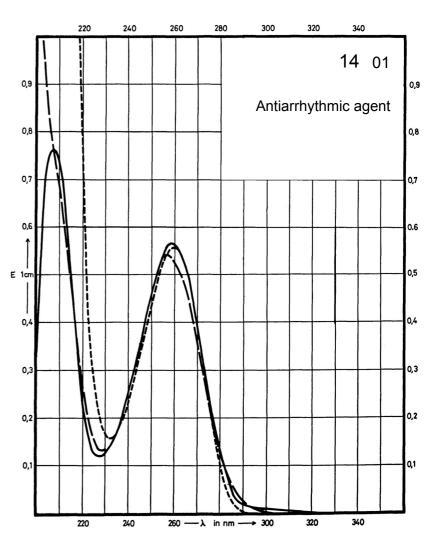

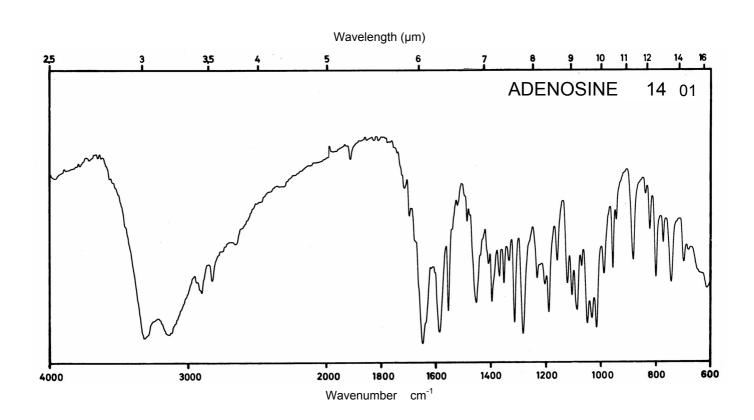

Name AZOSEMIDE


 M_{r}


370.8

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	339 nm 279 nm 237 nm	324 nm 282 nm 238 nm	Decom- position observed	328 nm 274 nm 236 nm
E ^{1%} 1cm	178 570 1453	162 488 1155		142 462 1187
ε	6600 21100 53900	6000 18100 42800		5300 17100 44000

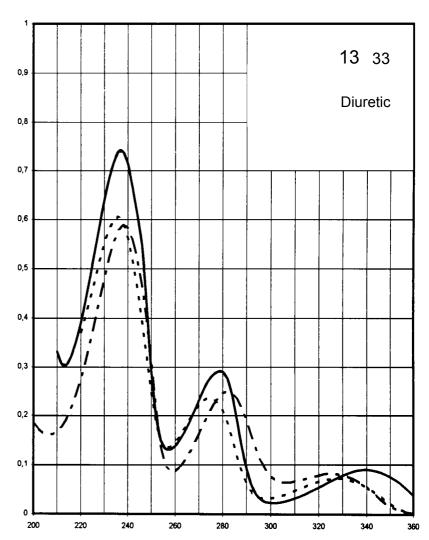

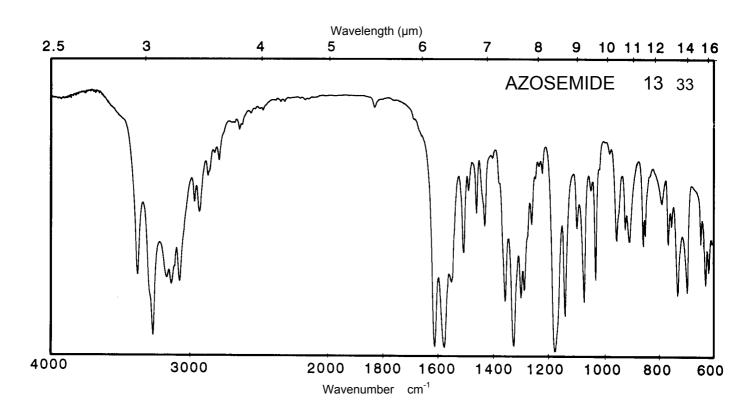



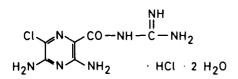

 M_r


302.1


Concentration

0.9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	361 nm 286 nm		361 nm 286 nm	366 nm 284 nm
E ^{1%} *	738 660		705 630	725 635
ε	19640 17560		18760 16760	19290 16900

^{*} Calculated on dried substance

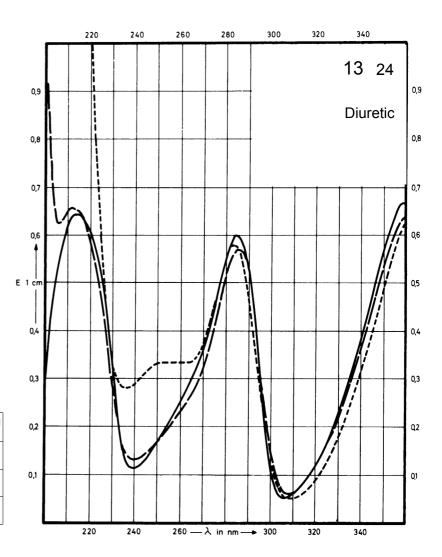

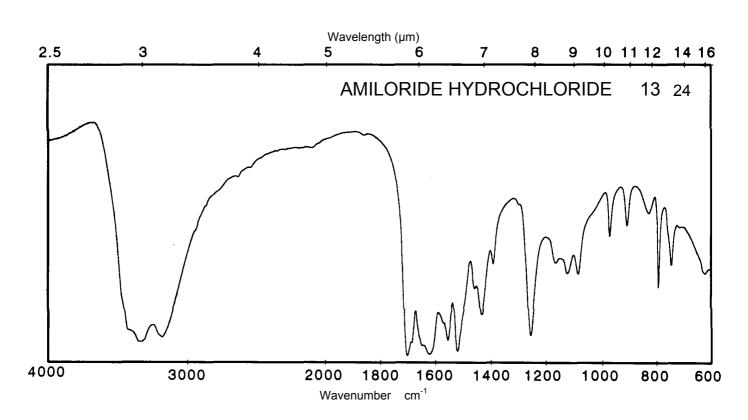



Name ACETAZOLAMIDE

 M_{r}


222.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm		263 nm	291 nm
E ^{1%} 1cm	433		450	551
ε	9630		10005	12250

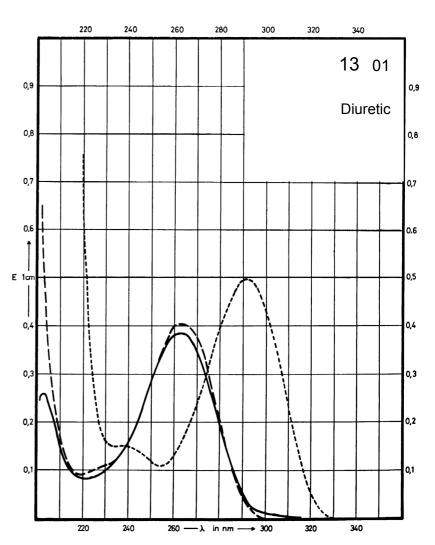

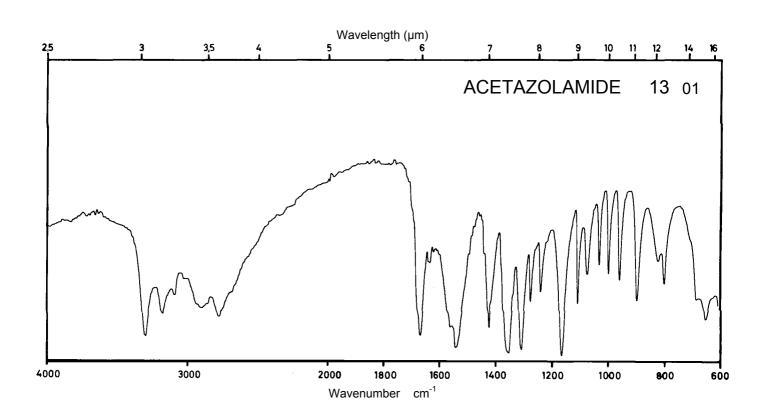



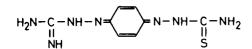


 M_r


237.3

Concentration

0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	440 nm		390 nm	460 nm
E ^{1%} 1cm	1040		1090	1100
ε	24680		25870	26100

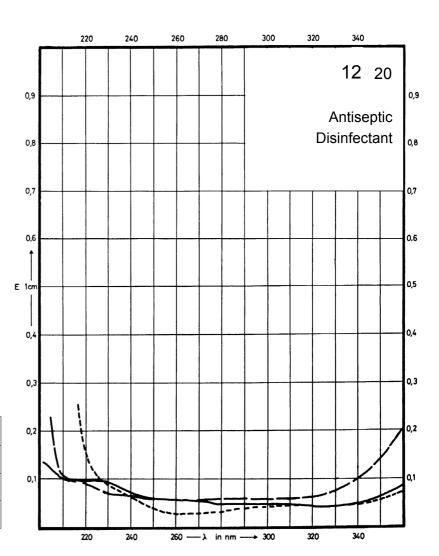


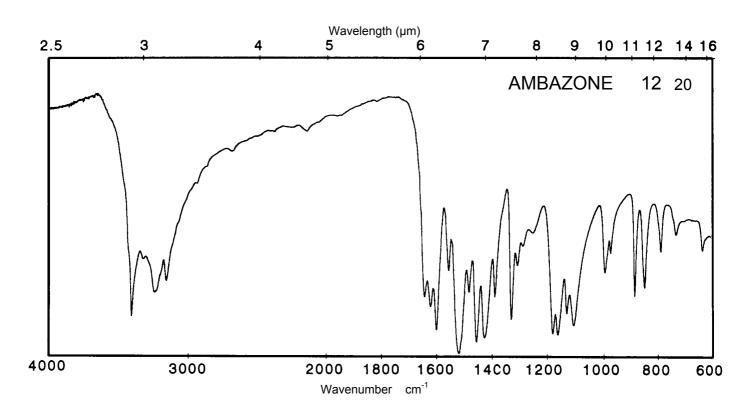


Name AZATHIOPRINE


 M_r 277.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm		280 nm	285 nm
E ^{1%} 1cm	633		621	560
ε	17600		17200	15500

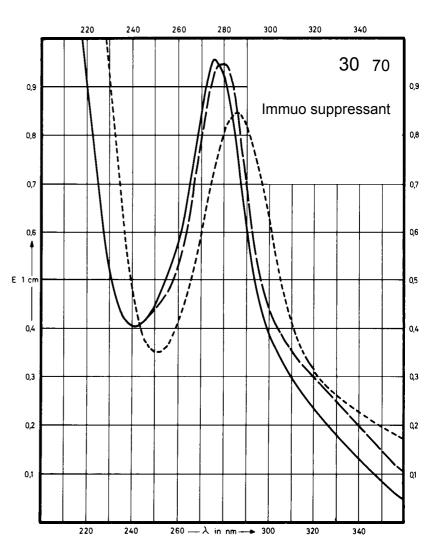

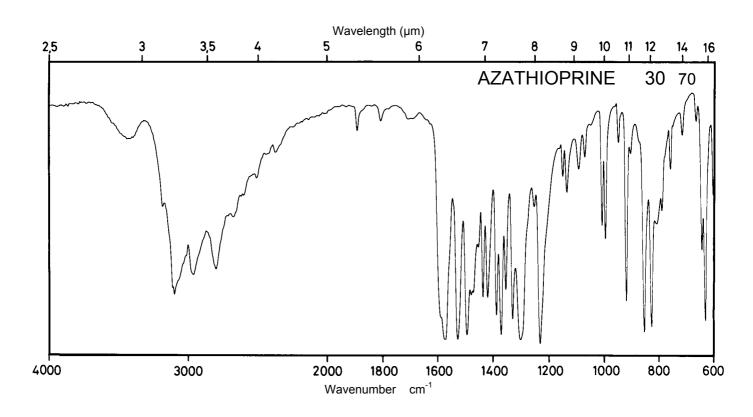



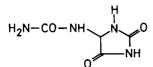



158.1

Concentration


1.2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				224 nm
E ^{1%} 1cm				375
ε				5900

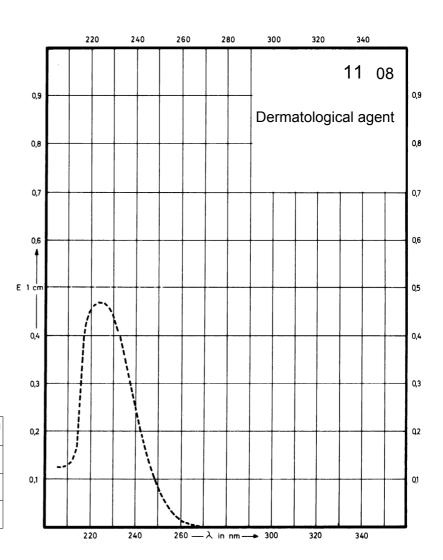

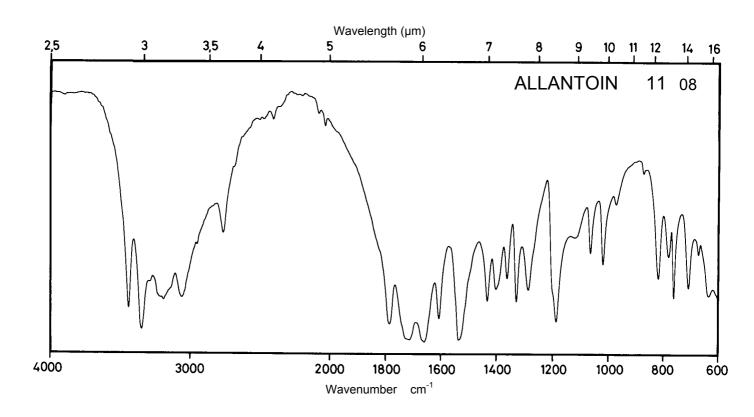



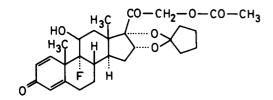



502.6

Concentration


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm			
E ^{1%} 1cm	308			
ε	15480			

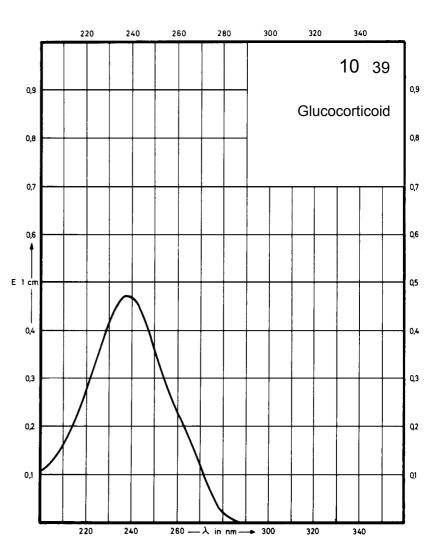

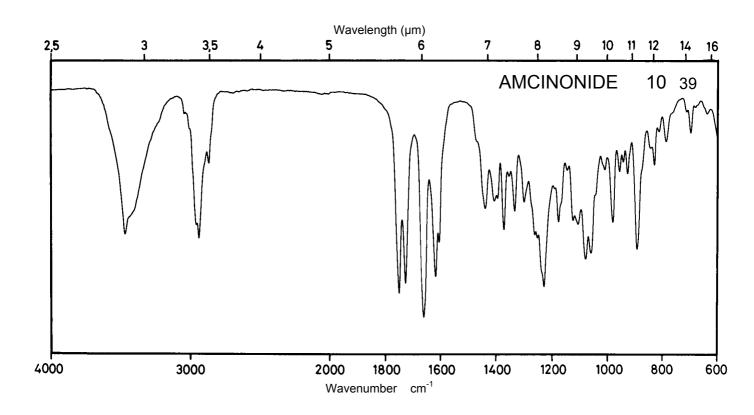


Name ALDOSTERONE

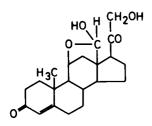


 M_r 360.5


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm	247 nm		
E ^{1%} 1cm	437	443		
3	15770	15980		

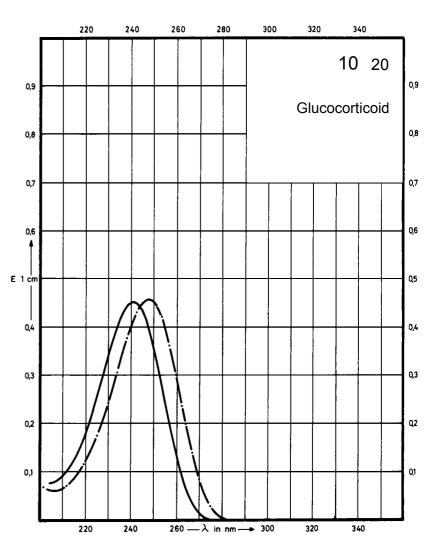

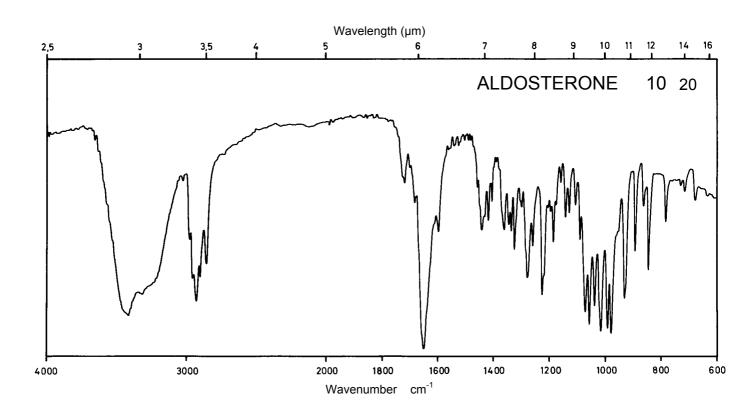



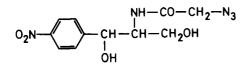



 M_r

295.3


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		278 nm	279 nm
E ^{1%} 1cm	339		326	303
ε	10000		9600	9000

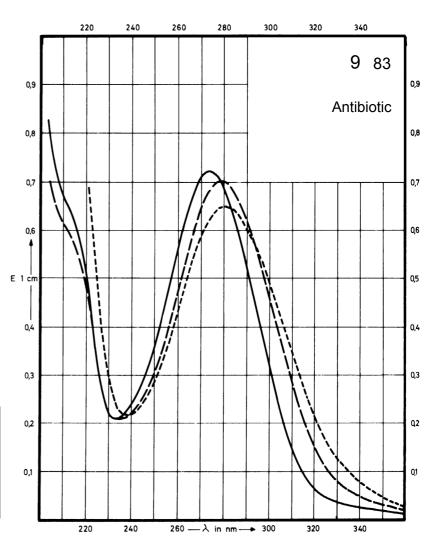

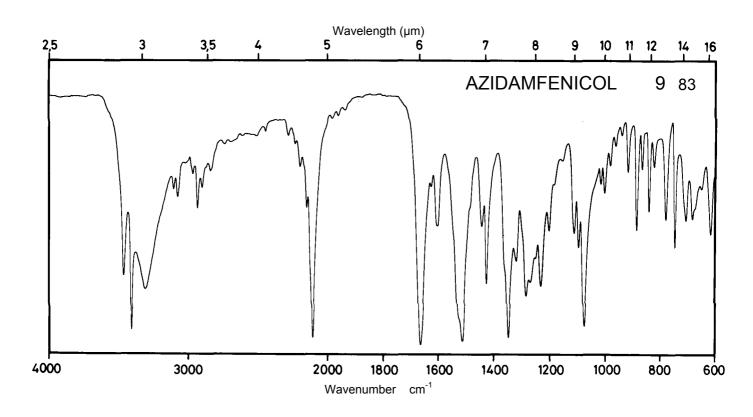



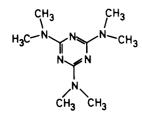



 $\,M_r\,$

210.3


Concentration 0.25 mg / 100 ml 0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	226 nm		240 nm	
E 1% 1cm	2400		1300	
ε	50000		27000	

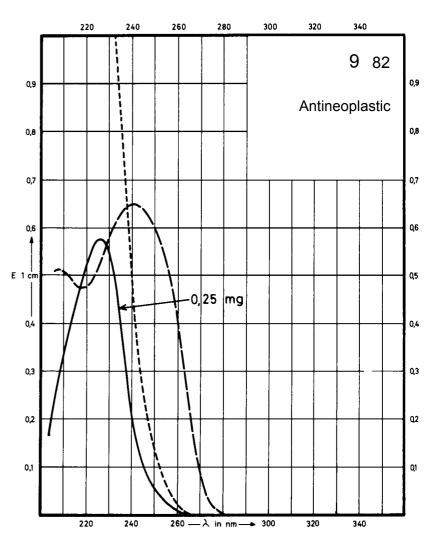

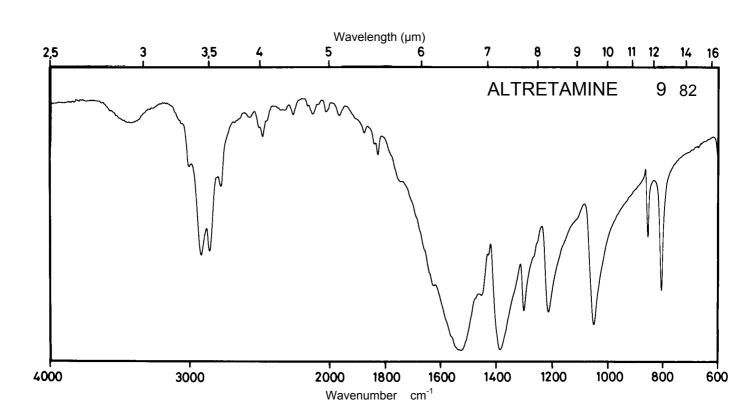


Name ACICLOVIR

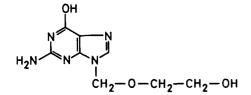


 M_r 225.2

Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	254 nm		255 nm	257 nm 264 nm
E ^{1%} 1cm	611		530	475 475
ε	13800		11900	10700 10700

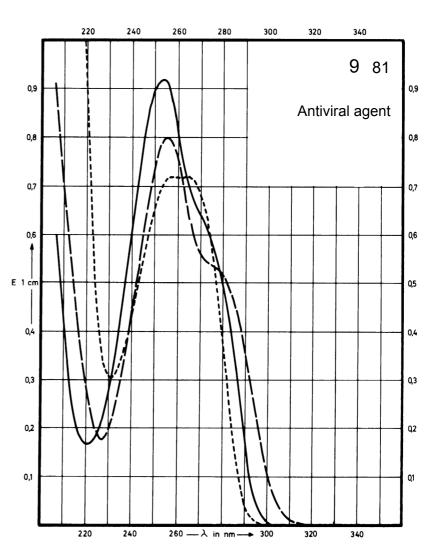

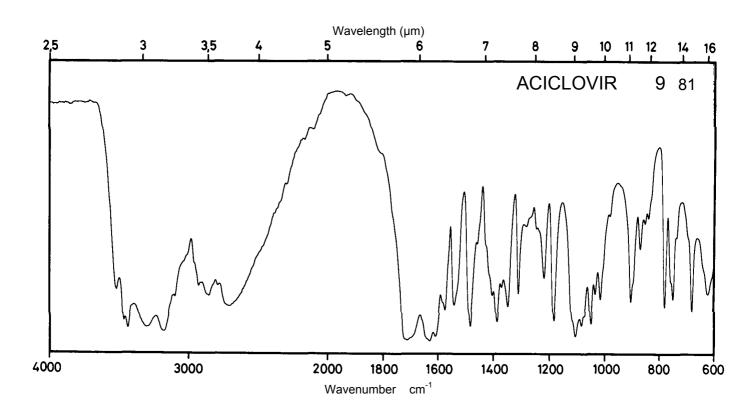

Name AMINOACRIDINE HYDROCHLORIDE

 M_{r}


248.7


Concentration

0.2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm		260 nm	260 nm
E ^{1%} _{1cm}	3500		3650	3260
ε	87050		90780	81080

Further Maxima (E 1%, 1cm)
Methanol: 422 nm (420); 399 nm (502); 380 nm (313)
0.1 M HCl: 422 nm (355); 399 nm (448); 380 nm (291)
0.1 M NaOH: 404 nm (350); 387 nm (307)

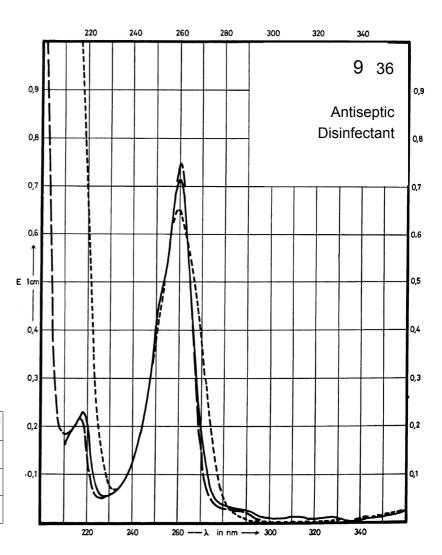

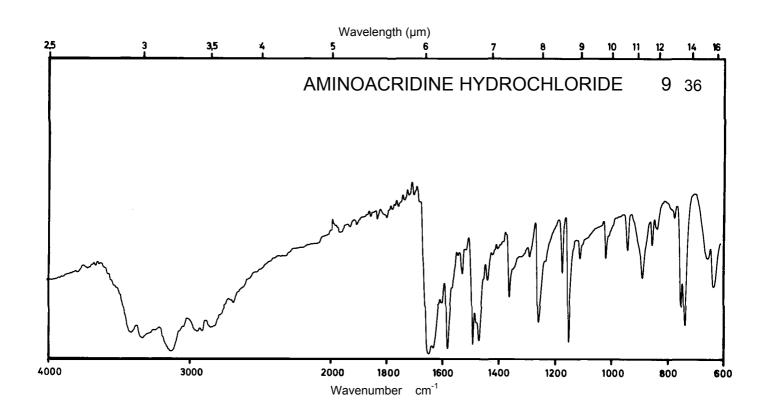

Name ARSANILIC ACID


O = AS OH


 M_{r}


217.1

Concentration 1 m


1 mg / 100 ml

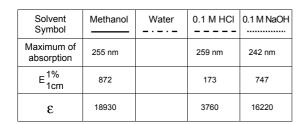

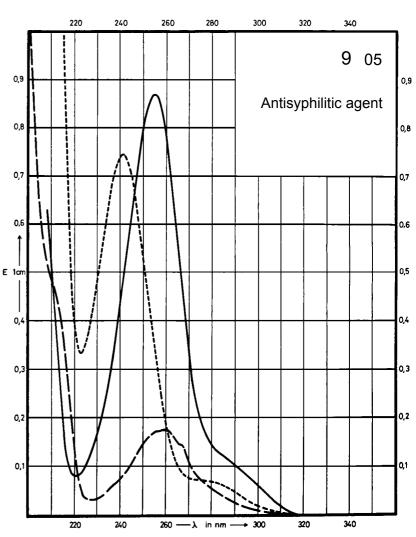

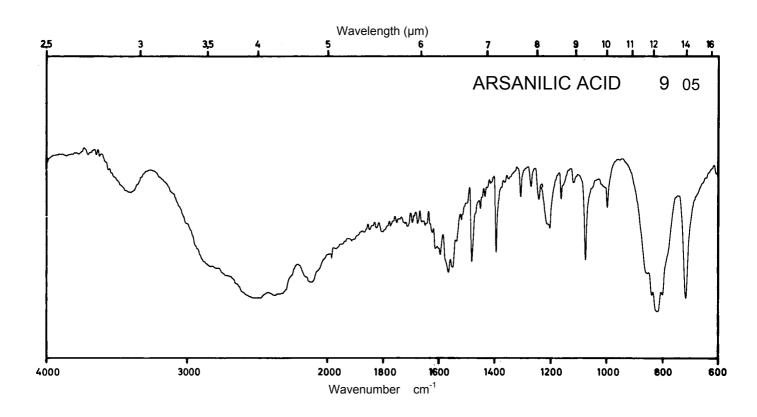

Name AMINOSALICYLIC ACID SODIUM SALT

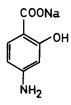


 M_{r}


175.1

Concentration


0.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	300 nm 266 nm		300 nm 233 nm	299 nm 264 nm
E ^{1%} 1cm	457 628		243 371	411 620
ε	8000 11000		4250 6500	7200 10860

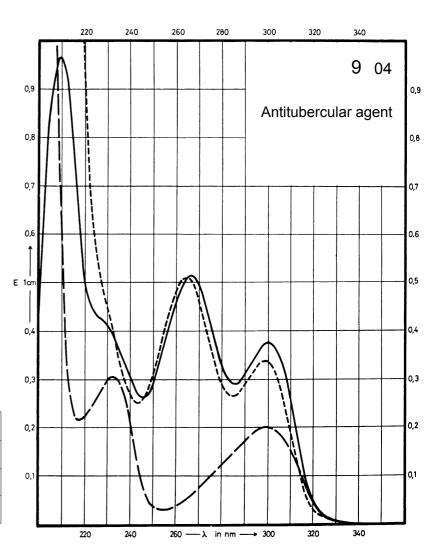

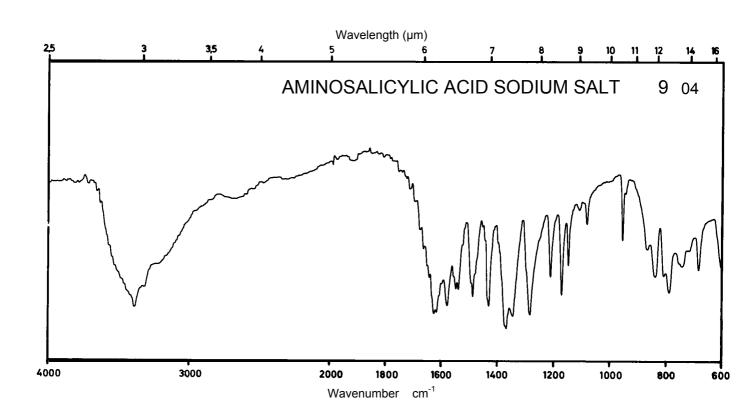




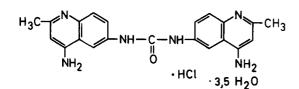


 M_r 472.0


Concentration 0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	345 nm 272 nm		342 nm 266 nm	
E ^{1%} 1cm	305 1603		253 1331	
ε	14400 75660		11940 62820	

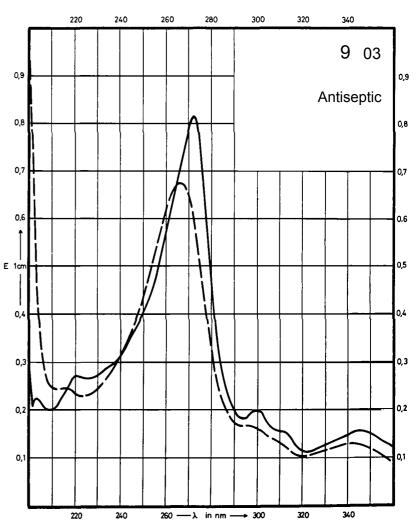

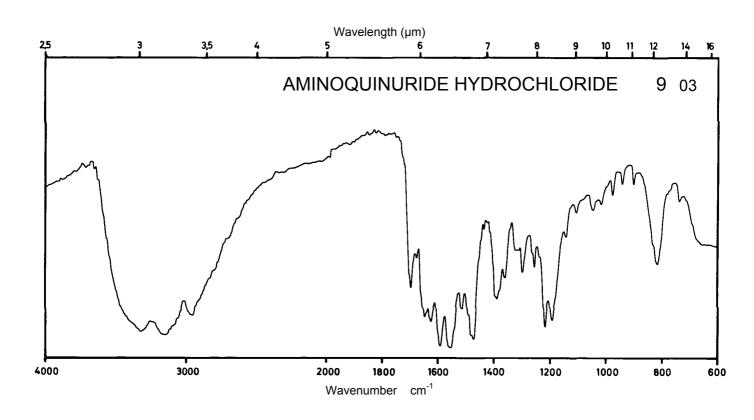




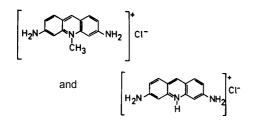



259.7 / 245.7

Concentration


0.25 mg / 100 ml

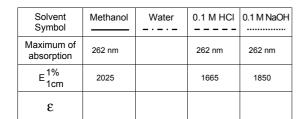

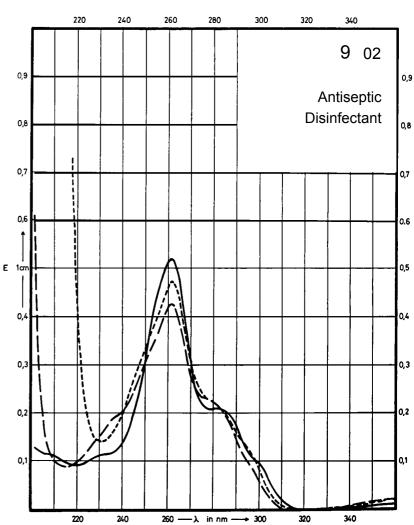

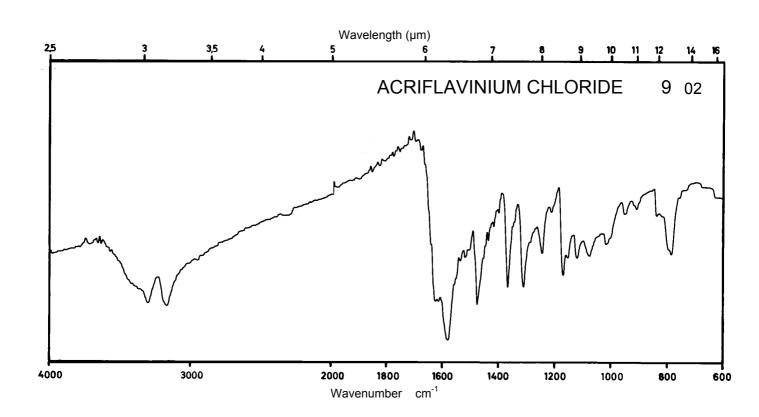

Name ACETARSOL

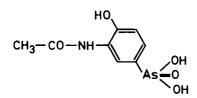



 M_{r}

275.1


Concentration 1.1 mg / 100 ml 5.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 249 nm		281 nm	299 nm 250 nm
E ^{1%} 1cm	150 399		94	177 567
ε	4130 10980		2590	4870 15600

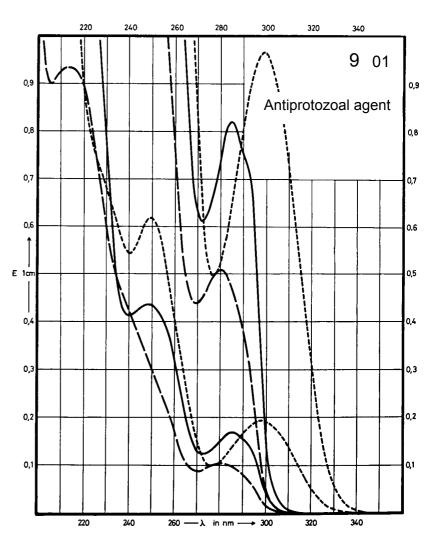

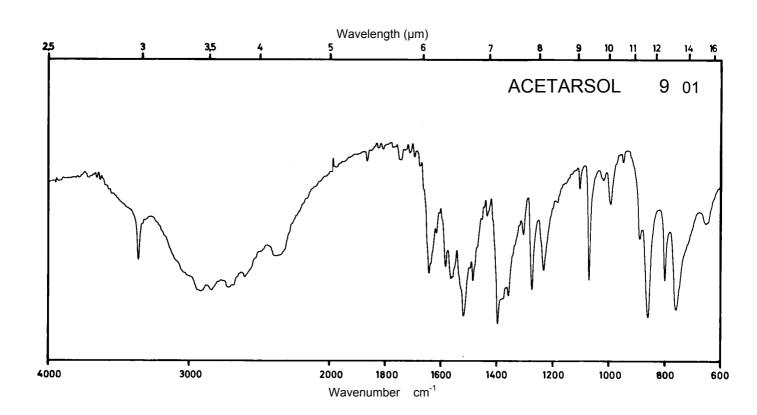



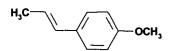



148.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm	257 nm	257 nm	257 nm
E ^{1%} 1cm	1556	1443	1424	1427
3	23100	21400	21100	21100

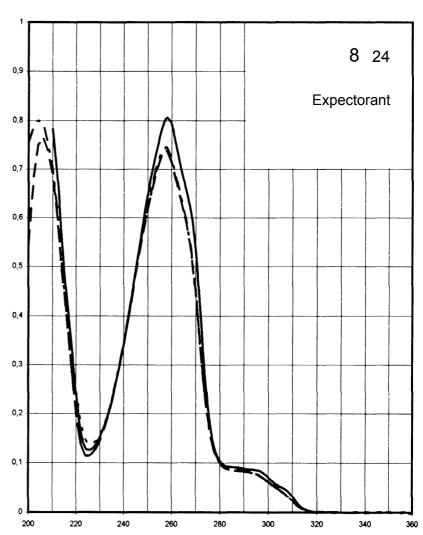

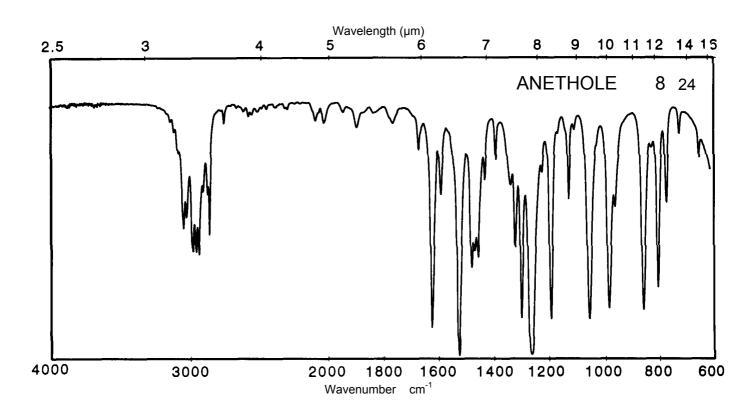

Name AMBROXOL HYDROCHLORIDE

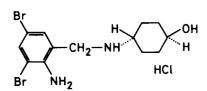


 M_{r}


414.6

Concentration


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	313 nm 248 nm		307 nm 244 nm	299 nm 241 nm
E ^{1%} _{1cm}	79 282		72 241	64 240
ε	3300 11700		3000 10000	2700 9900

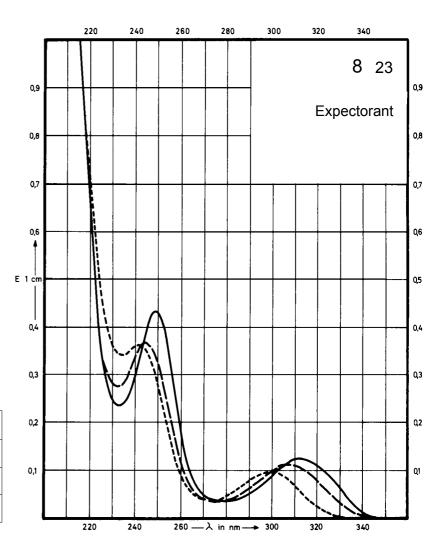

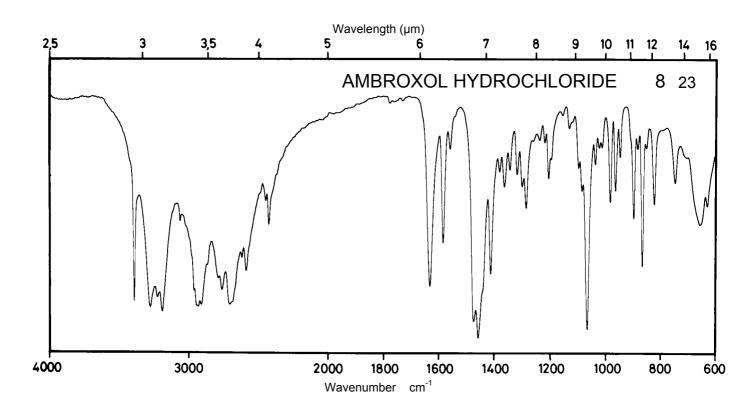




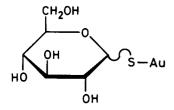


392.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

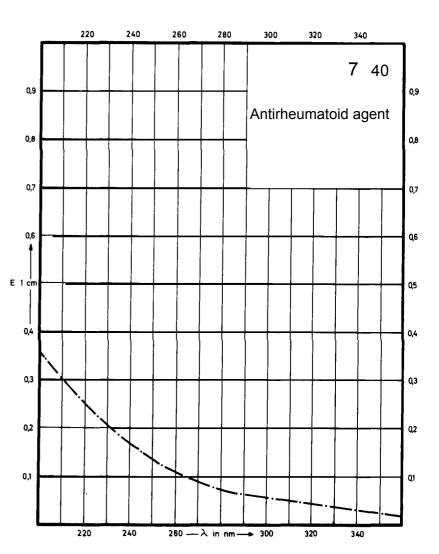

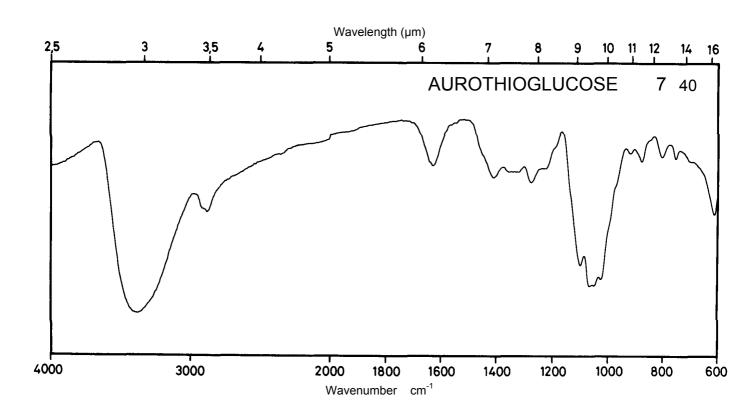




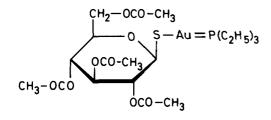


678.5


Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	221 nm			
E ^{1%} 1cm	86			
ε	5840			

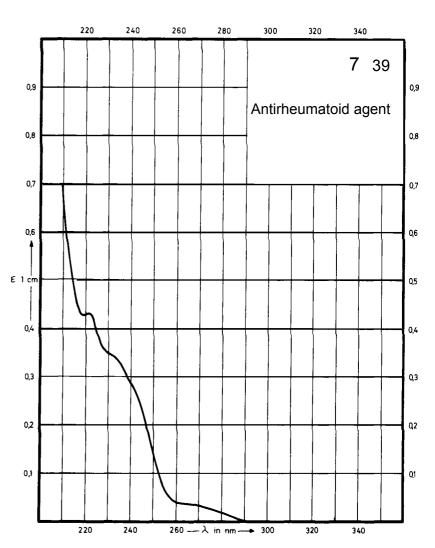

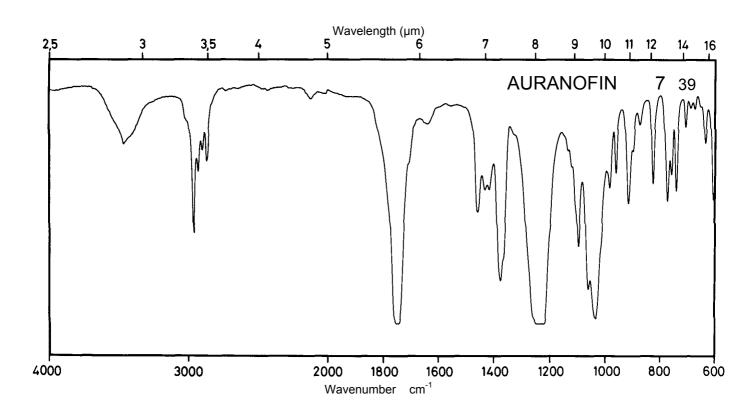




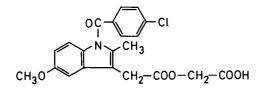


 M_r 415.8

Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	319 nm	319 nm 263 nm	319 nm 263 nm	281 nm
E ^{1%} 1cm	158	173 455	139 368	165
ε	6570	7190 18920	5780 15300	6870

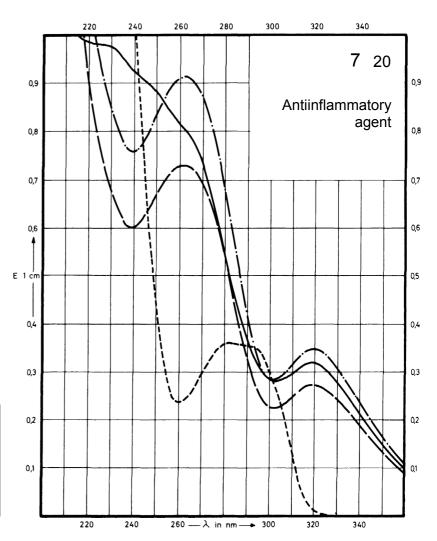

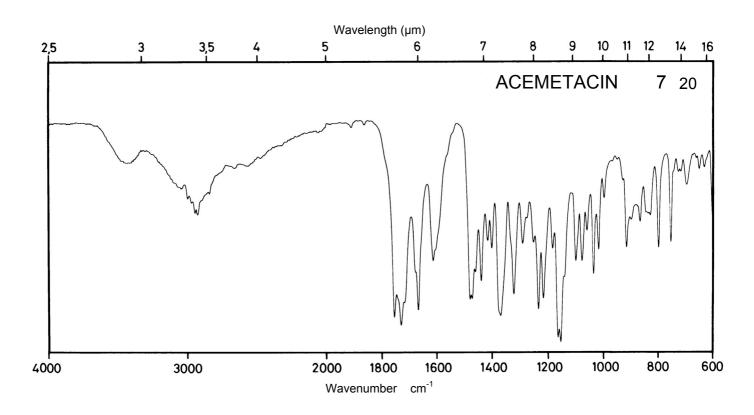

Name AZAPROPAZONE DIHYDRATE


 M_{r}


336.4

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	251 nm		251 nm	325 nm 256 nm
E ^{1%} 1cm	1300		1040	223 1080
ε	43730		34980	7500 36420

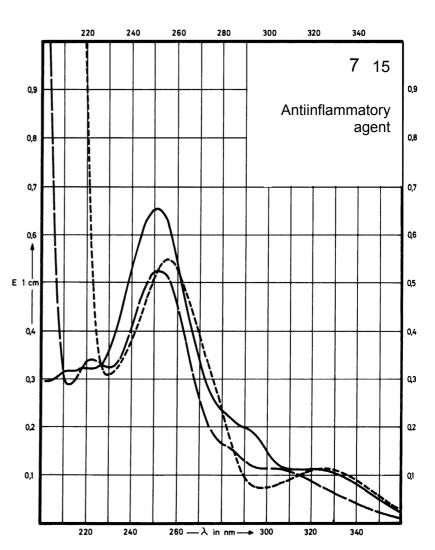

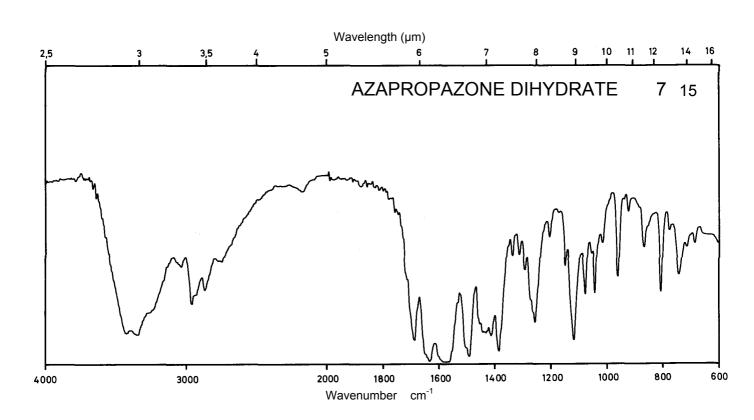




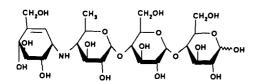


645.6


Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

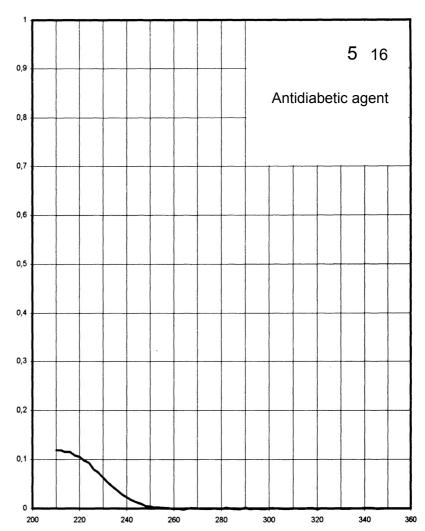

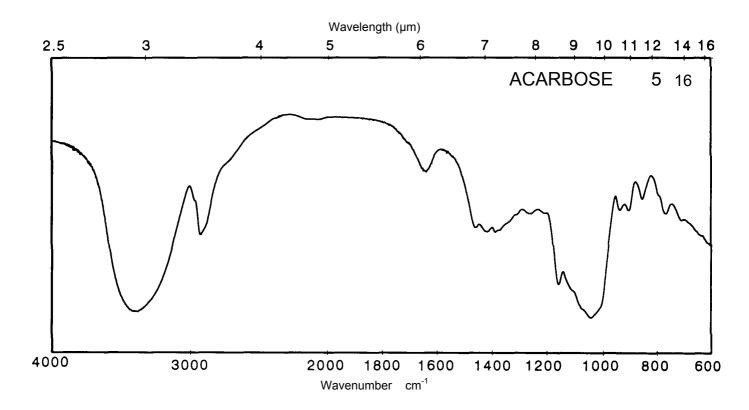




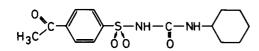


324.4


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	246 nm		246 nm	247 nm
E ^{1%} 1cm	424		508	413
3	13750		16480	13400

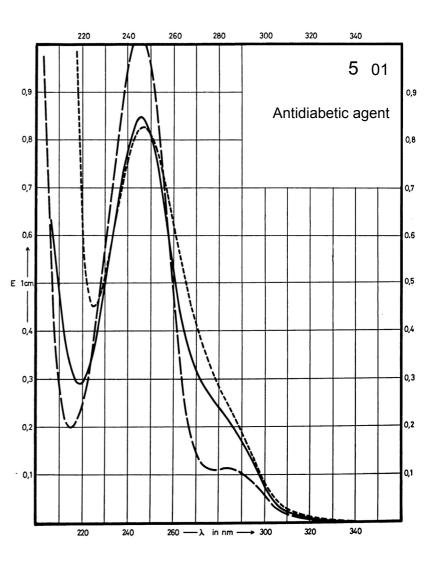

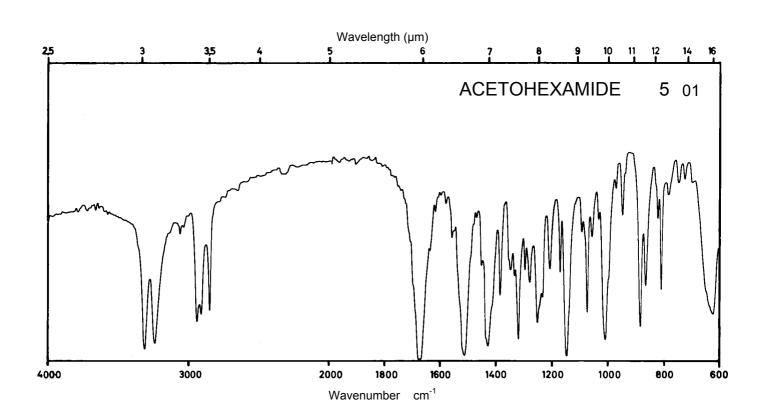



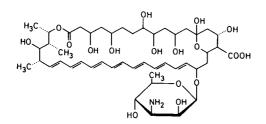



924.1

Concentration


1.06 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				406 nm 385 nm 366 nm
E 1% 1cm				1075 1070 690
ε				99300 98600 63600

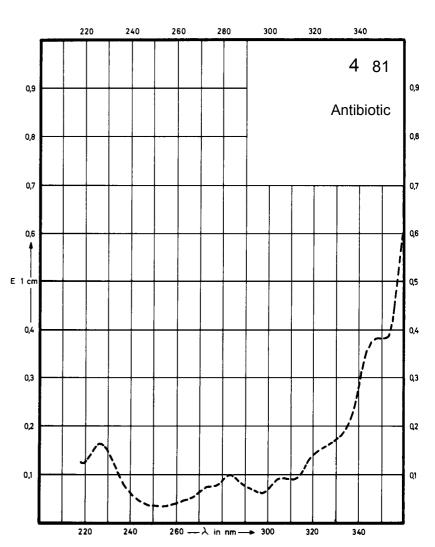

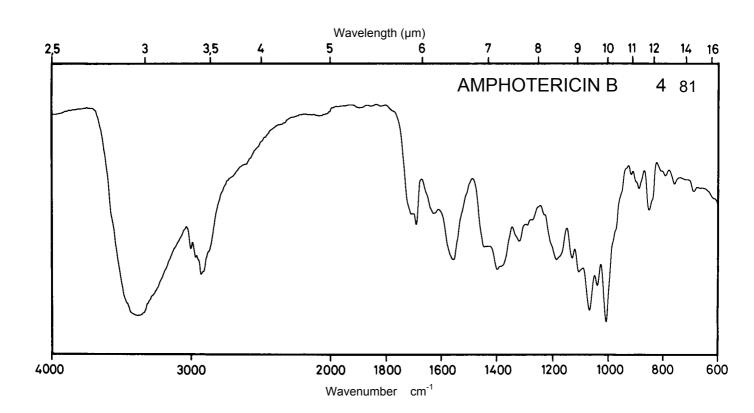




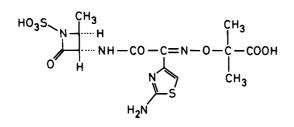


 M_r 435.4


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	235 nm	285 nm	282 nm 266 nm	
E ^{1%} 1cm		206	279 277	
3		8990	12170 12070	

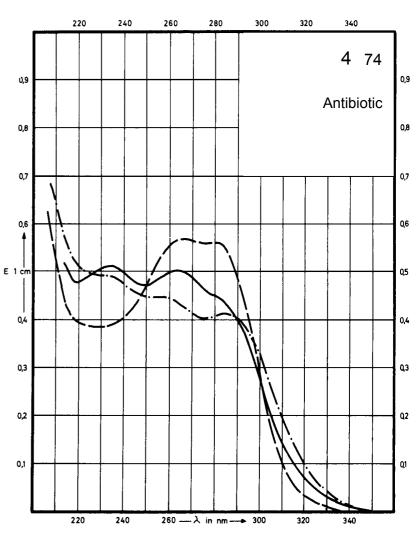

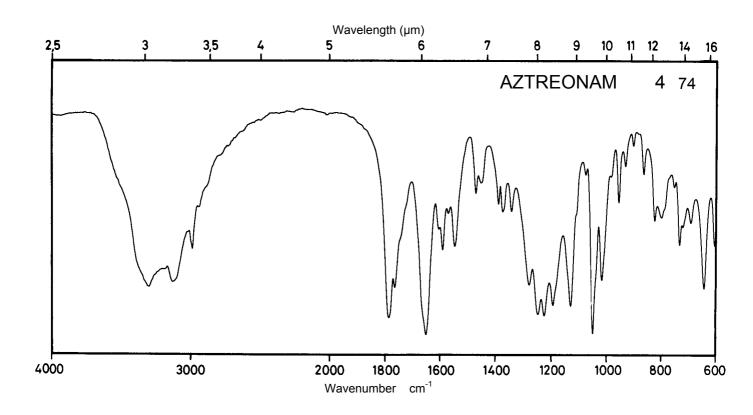




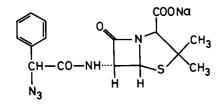


397.4


Concentration 60 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm	256 nm		
E ^{1%} 1cm	7.9	8.6		
3	312	342		

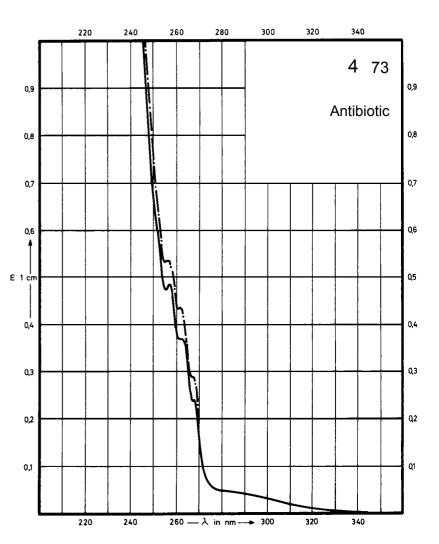

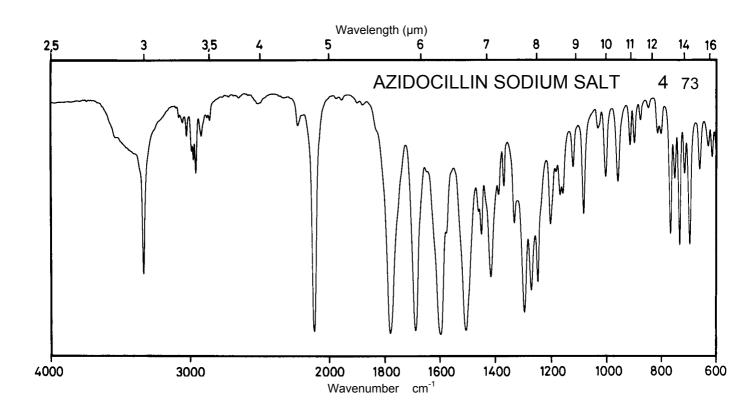



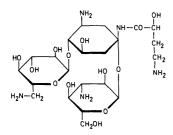



585.6

Concentration 5


50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

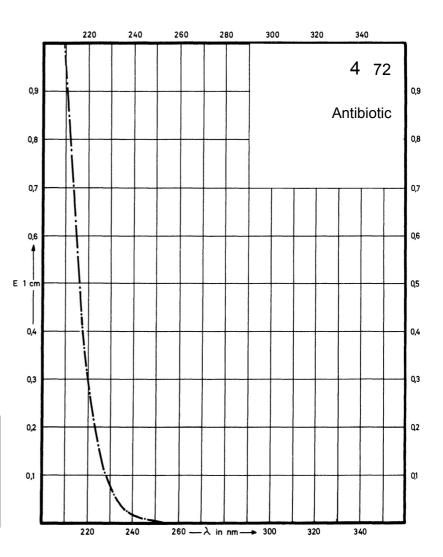

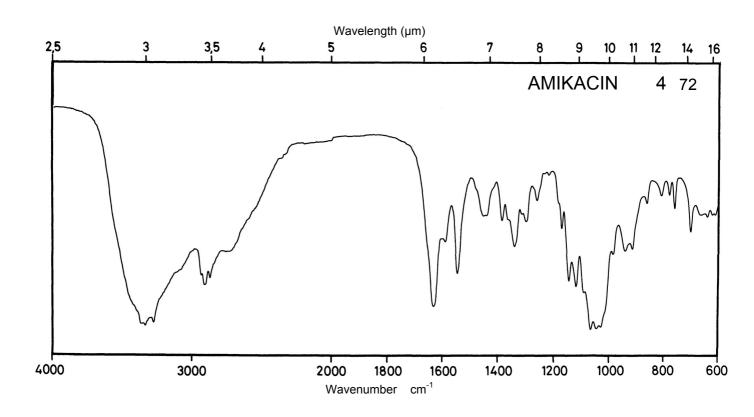


Name APALCILLIN SODIUM SALT

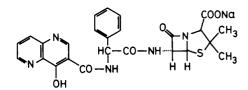



543.5


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 255 nm		335 nm 264 nm	317 nm 264 nm
E ^{1%} 1cm	75 261		36 310	66 365
ε	4100 14200		2000 16900	3600 19800

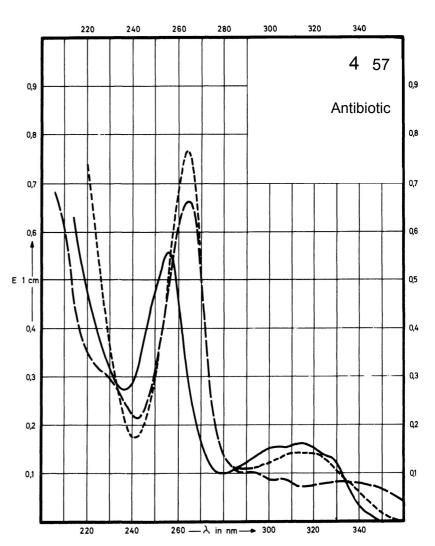

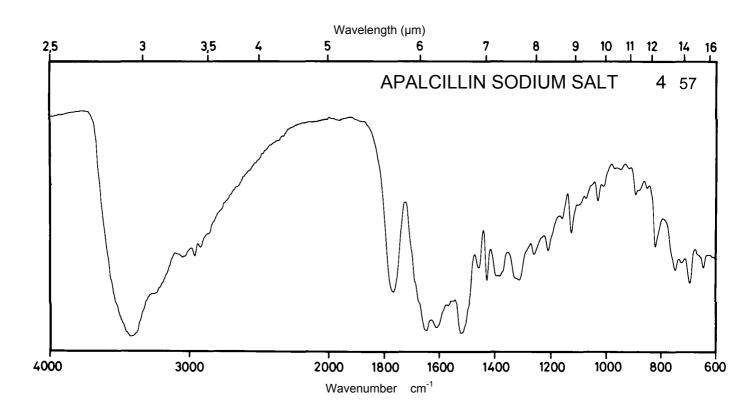

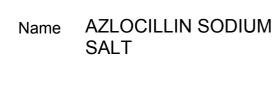


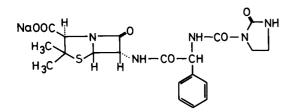


483.5


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 258 nm	257 nm		
E ^{1%} 1cm	7.4 10.8	13.1		
ε	360 520	634		

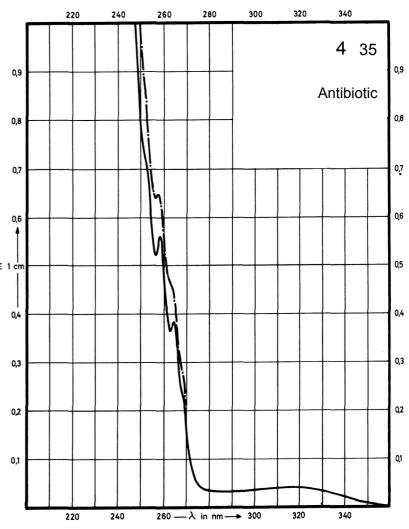

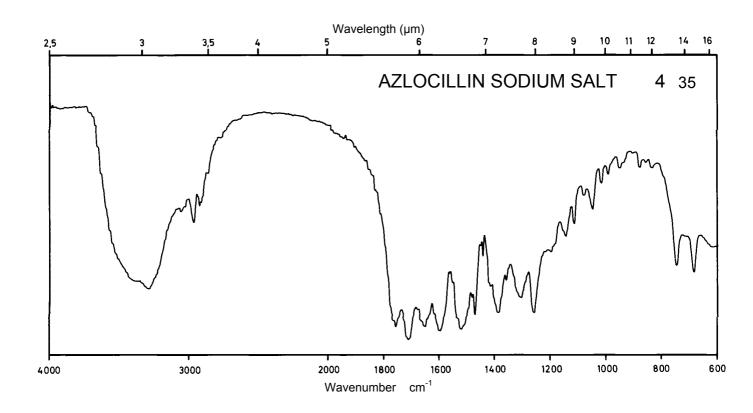




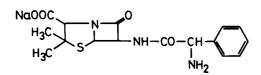


371.4


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 259 nm		269 nm 262 nm 257 nm	
E ^{1%} 1cm	4.6 7.0		5.6 7.9 8.5	
ε	170 260		208 295 316	

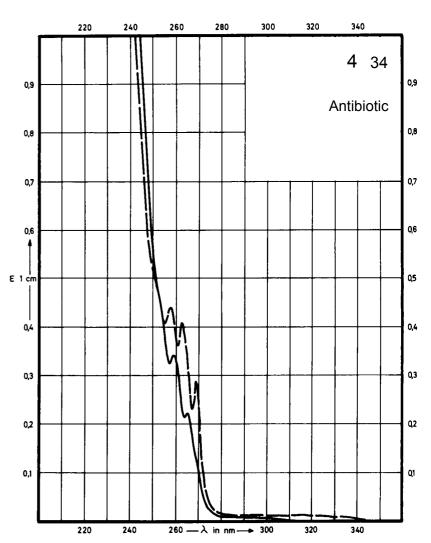

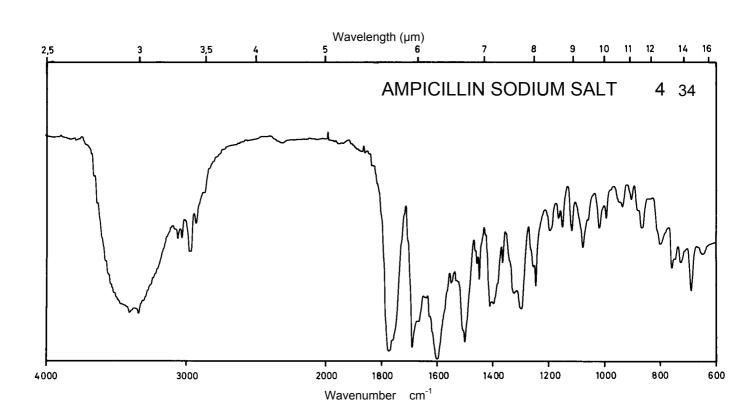



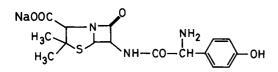



 M_{r}

387.4


Concentration 3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 229 nm		272 nm 228 nm	
E ^{1%} 1cm	36 236		27 232	
3	1390 9160		1050 9000	

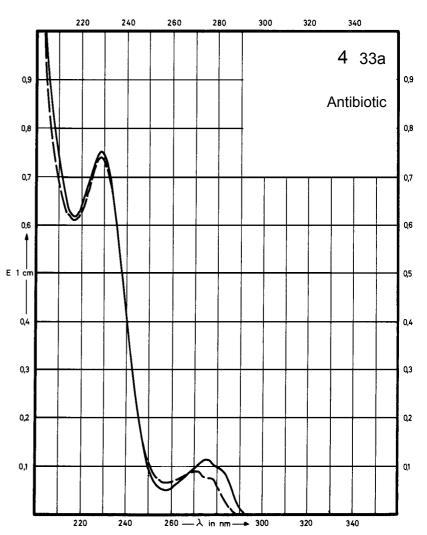

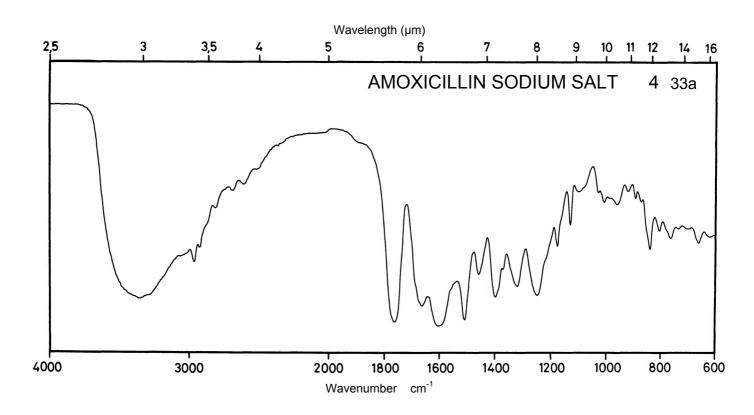




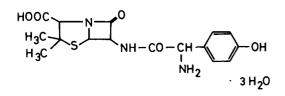


 M_r 419.5


Concentration 2.5 mg / 100 ml 5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			273 nm 230 nm	291 nm 247 nm
E ^{1%} 1cm			26 225	64 286
ε			1070 9450	2670 11980

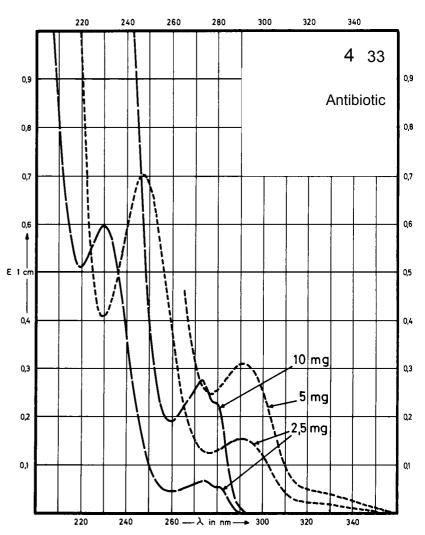

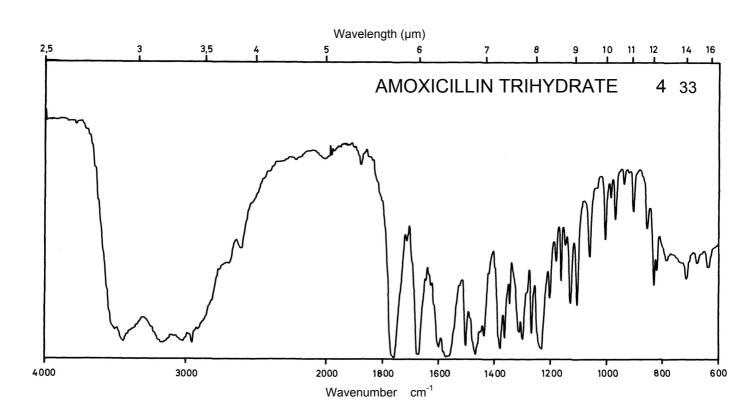




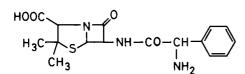


 M_r 349.4


Concentration 50 mg / 100 ml

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	267 nm	267 nm	267 nm	Decom-
	261 nm	261 nm	261 nm	position
absorption	256 nm	256 nm	256 nm	observed
1%	5.4	6.3	6.5	
E ^{1%} 1cm	7.6	8.9	9.0	
TCM	8.9	9.9	9.7	
	185	220	225	
3	265	310	315	
	240	250	240	1

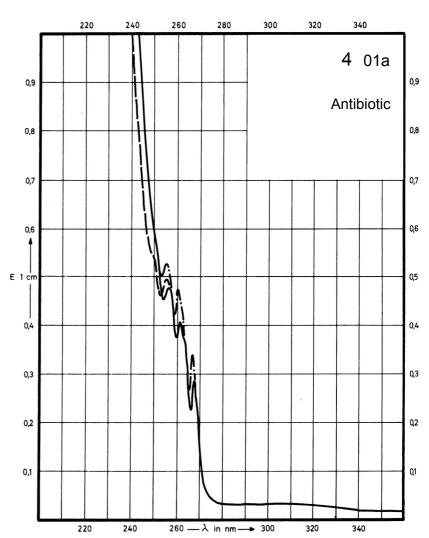

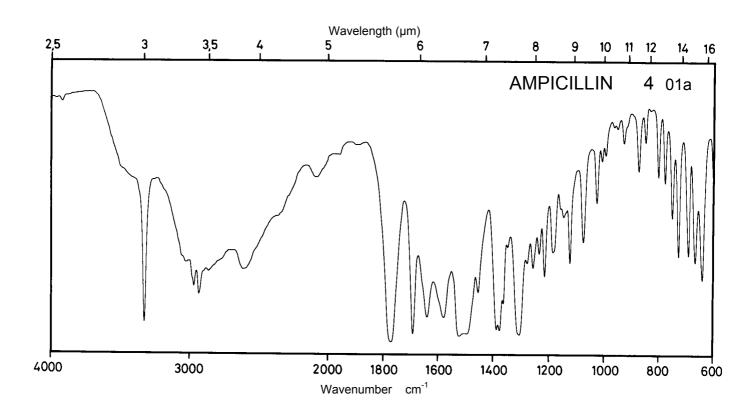



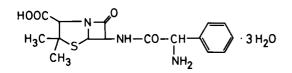



 M_r

Concentration 60 mg / 100 ml


403.5

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	267 nm	267 nm	267 nm	Decom-
	261 nm	261 nm	261 nm	position
absorption	256 nm	256 nm	256 nm	observed
_1%	4.6	5.5	5.6	
F	6.6	7.7	7.8	
-1cm	7.7	8.6	8.4	
	185	220	225	
3	265	310	315	
	310	350	340	

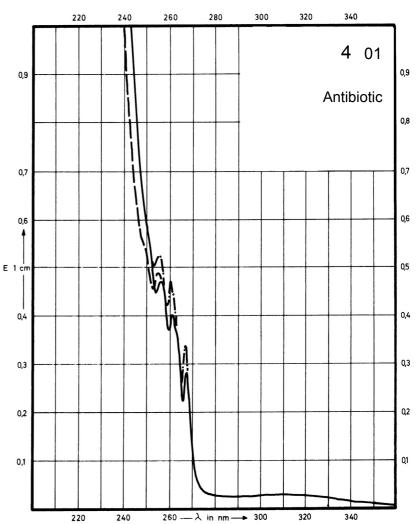

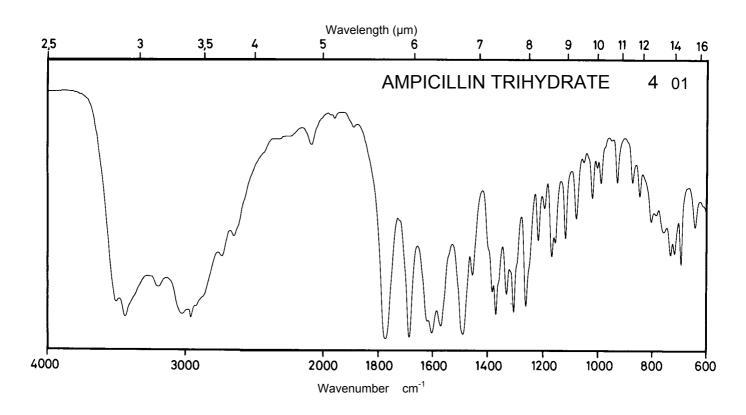




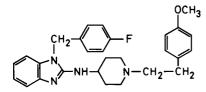


 M_r 458.6


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 250 nm		282 nm 276 nm	
E ^{1%} 1cm	268 200		272 289	
ε	12300 9150		12500 13300	

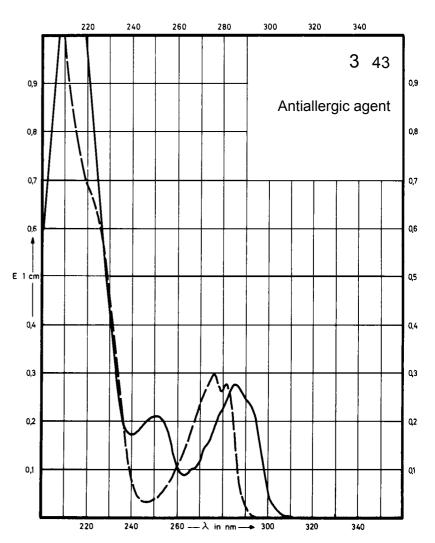

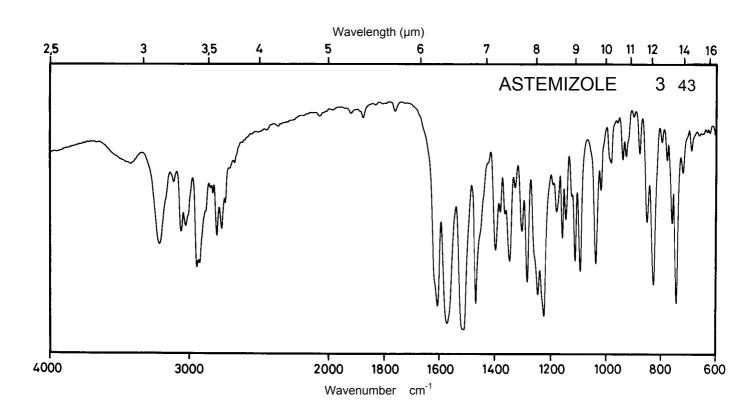



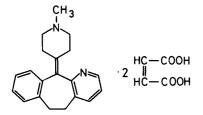



522.6

Concentration


2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm		283 nm	
E ^{1%} 1cm	127		177	
3	6660		9200	

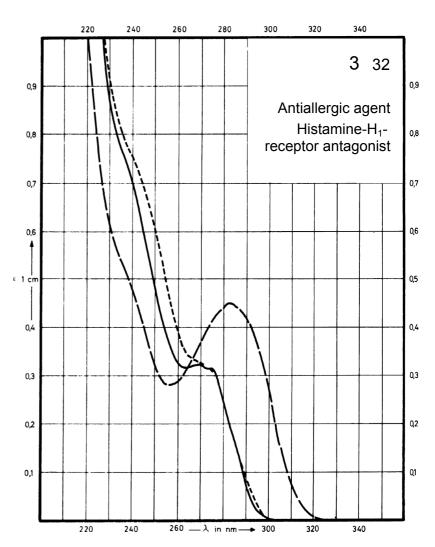

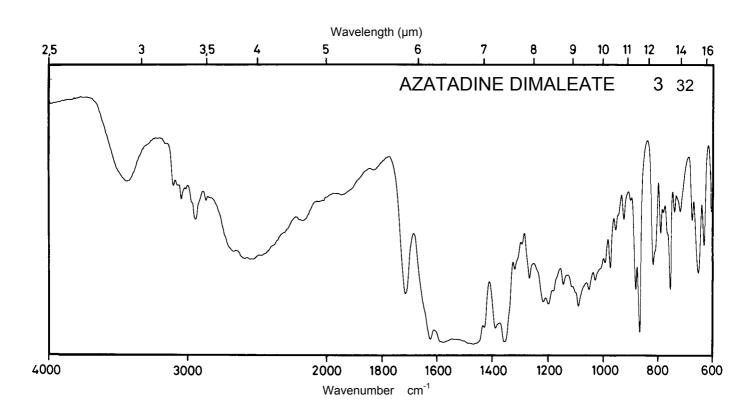



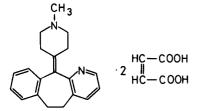



 M_r

Concentration 2.5 mg / 100 ml


522.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm		283 nm	
E ^{1%} 1cm	127		177	
3	6660		9200	

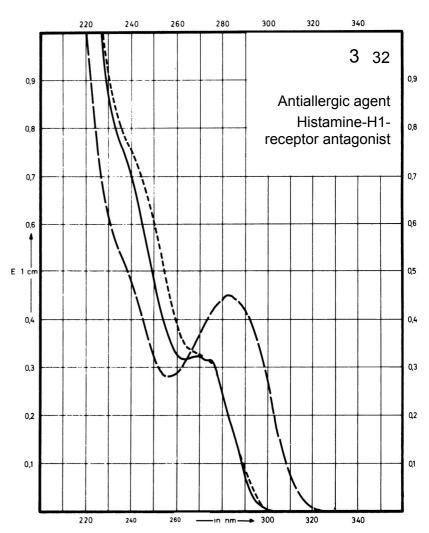

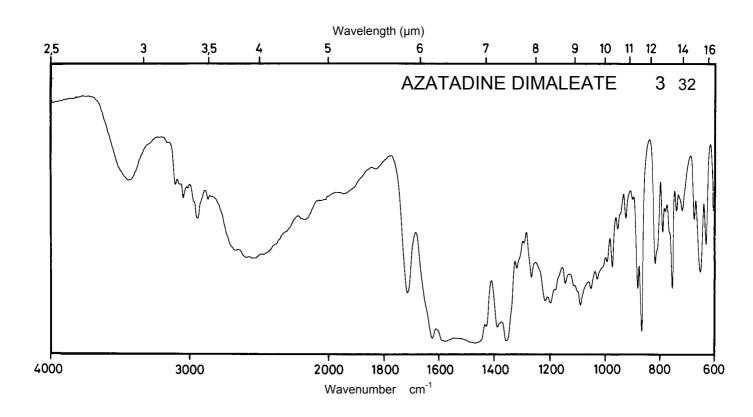



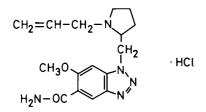



 M_r

351.8


Concentration 0.4 mg / 100 ml 3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm 222 nm	295 nm	295 nm	286 nm
E ^{1%} 1cm	178 1030	183	182	277
ε	6260 36240	6440	6400	9740

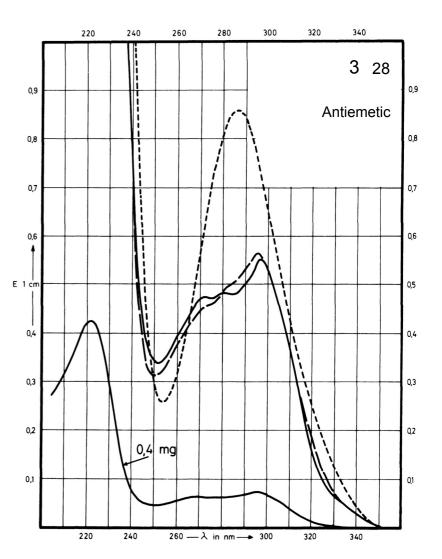

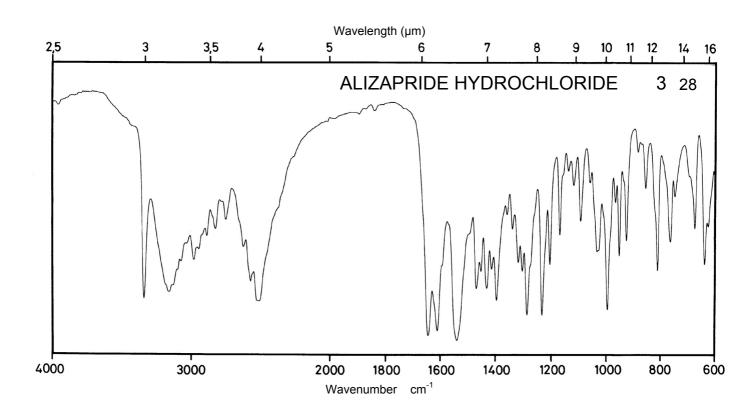




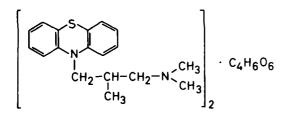


746.9


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	305 nm 254 nm		300 nm 251 nm	
E ^{1%} 1cm	112 840		99 760	
ε	8370 62740		7390 56760	

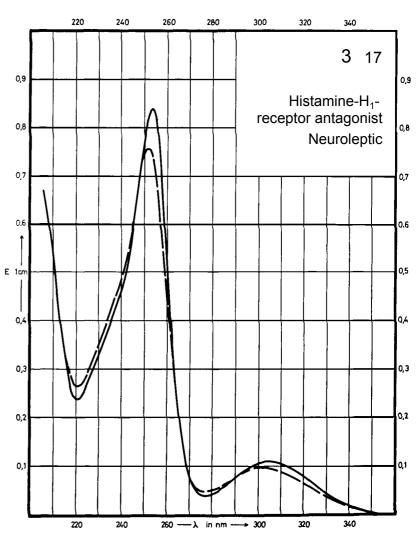

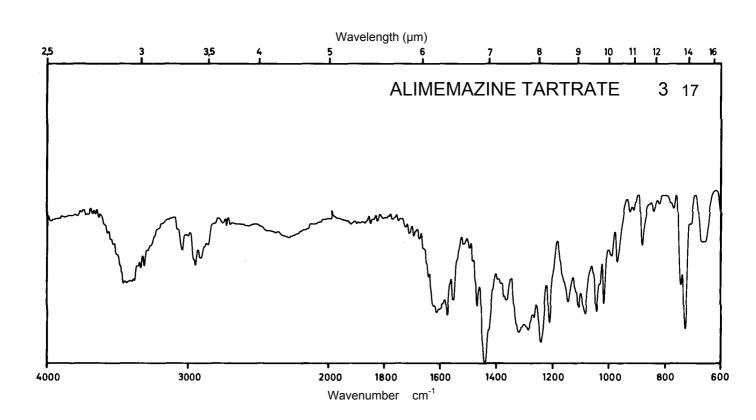

Name ANTAZOLINE HYDROCHLORIDE

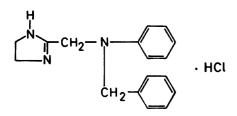



 M_r

Concentration 1 mg / 100 ml 10 mg / 100 ml


301.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	292 nm 240 nm		292 nm 241 nm	295 nm 248 nm
E ^{1%} 1cm	66 526		68 502	76 509
ε	1990 15870		2050 15150	2290 15360

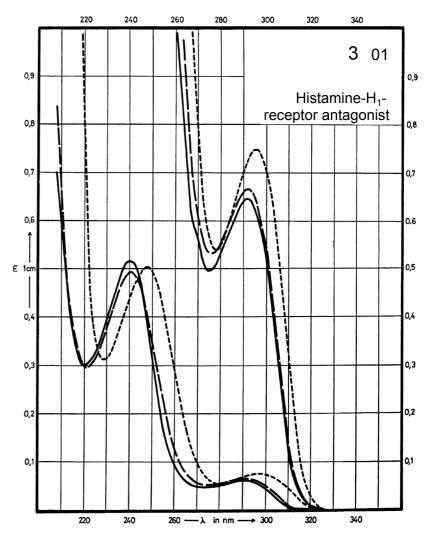

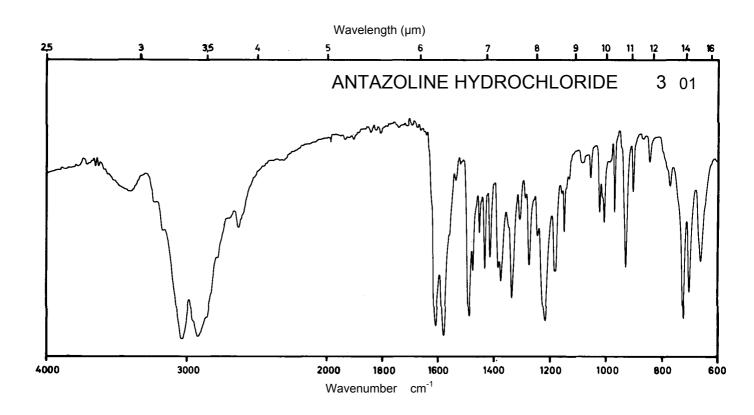





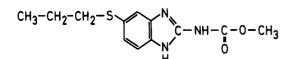


 M_{r} 265.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm		292 nm	309 nm
E ^{1%} 1cm	437		370	742
3	11580		9800	19700

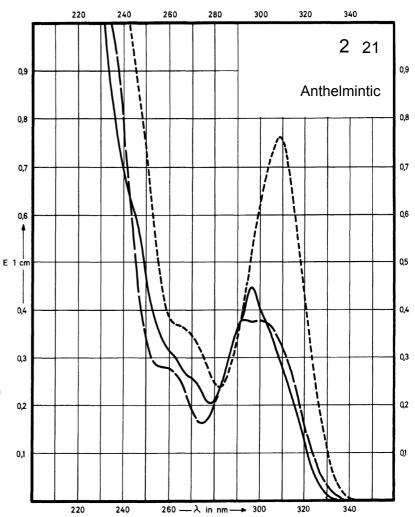

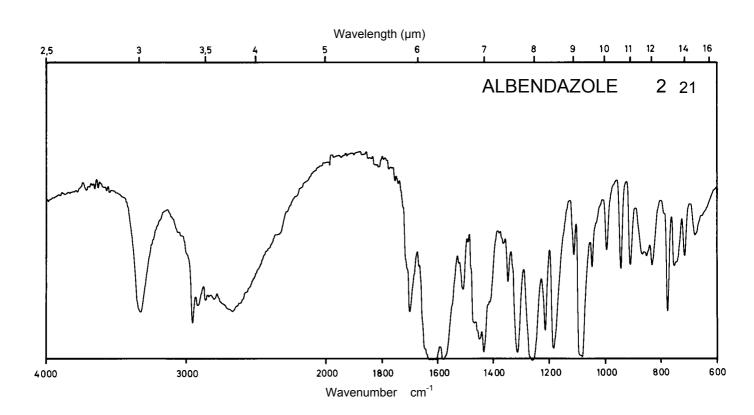




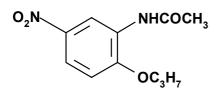


 M_r


238.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	297 nm 253 nm	319 nm 254 nm 227 nm	320 nm 254 nm 227 nm	322 nm 254 nm
E ^{1%} 1cm	323 633	353 419 423	353 421 425	342 421
ε	7690 15100	8400 10000 10000	8400 10000 10000	8150 10000

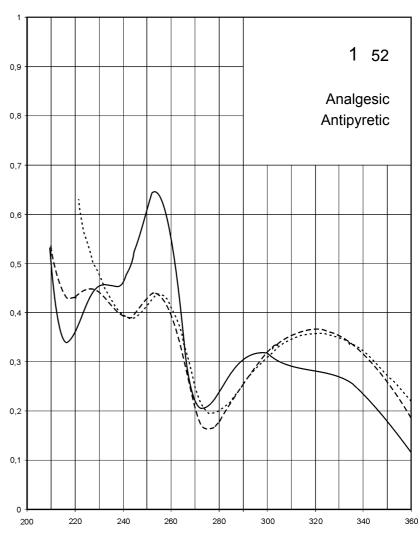

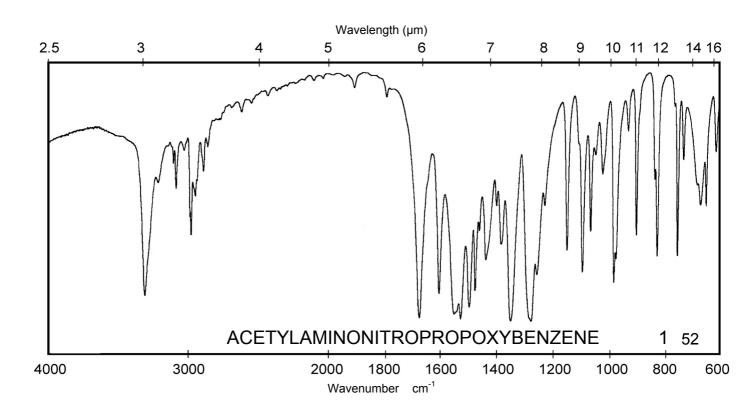




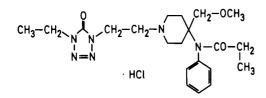



 M_r


453.0

Solvent Symbol	Propanol-2	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm			
E ^{1%} 1cm	5.8			
ε	263			

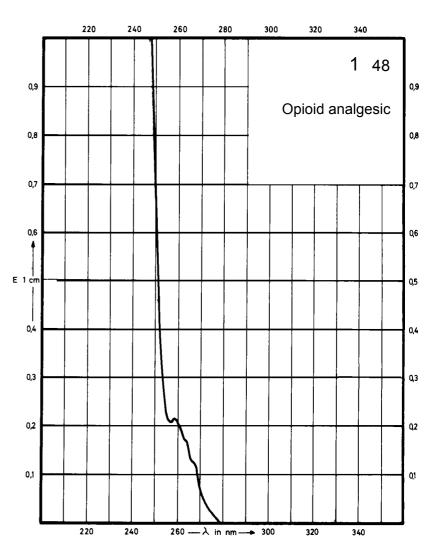

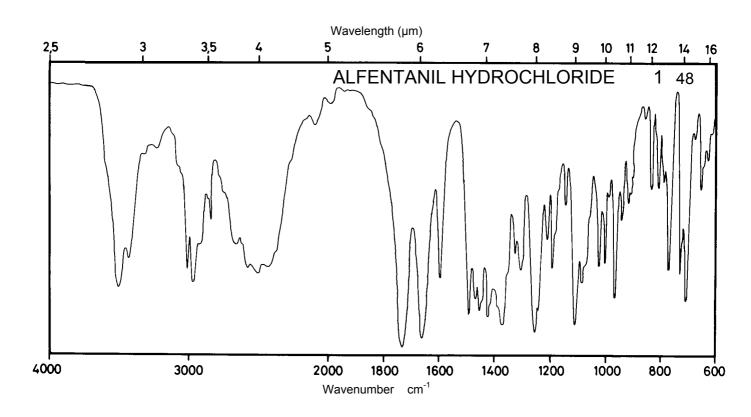

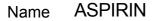


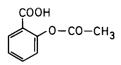


 M_r 180.2


Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 226 nm	269 nm	276 nm 227 nm	Decom- position observed
E ^{1%} 1cm	65 488	37	66 480	
ε	1170 8790	670	1190 8650	

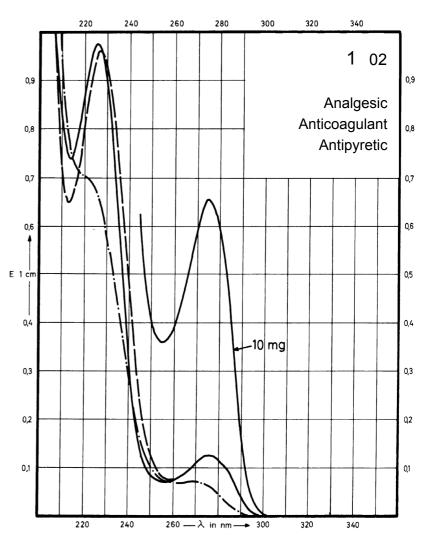

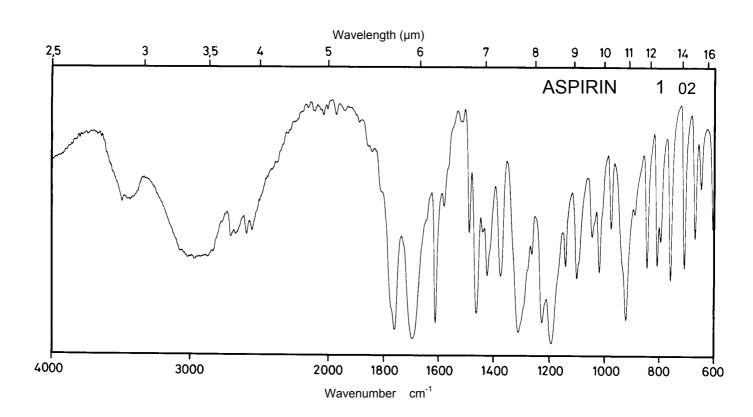



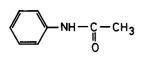



135.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm		237 nm	239 nm
E ^{1%} 1cm	1060		815	820
ε	14330		11020	11090

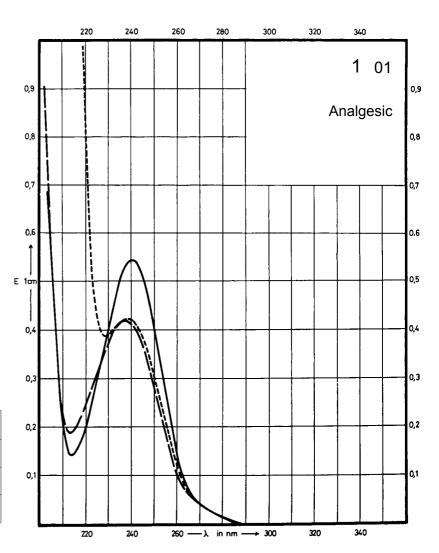

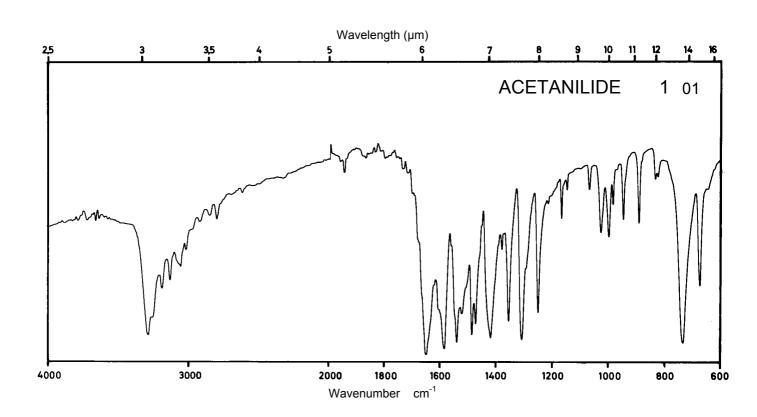



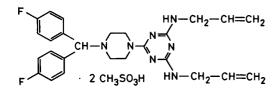



669.8


Concentration


0.76 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm 227 nm		228 nm	
E ^{1%} 1cm	522 500		630	
ε	35000 33500		42200	

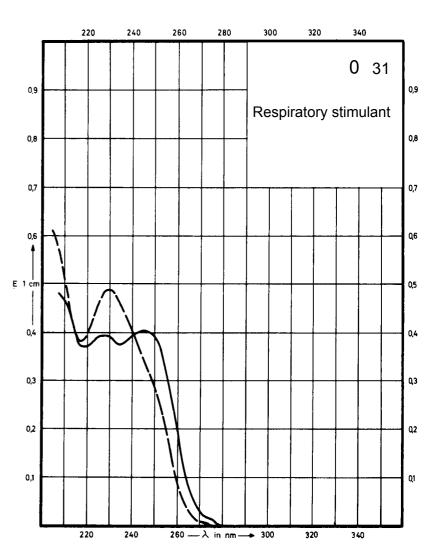

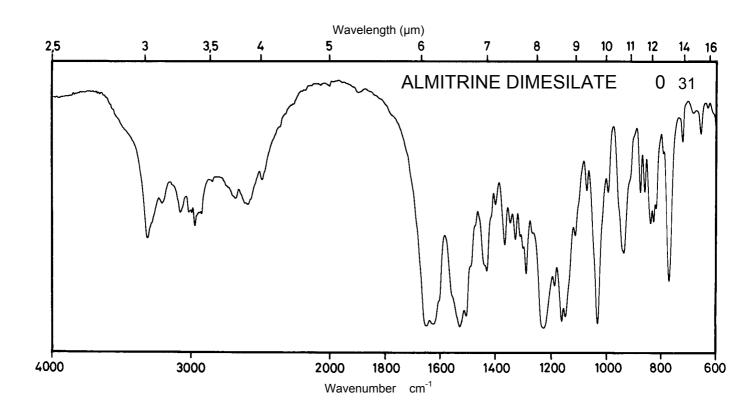

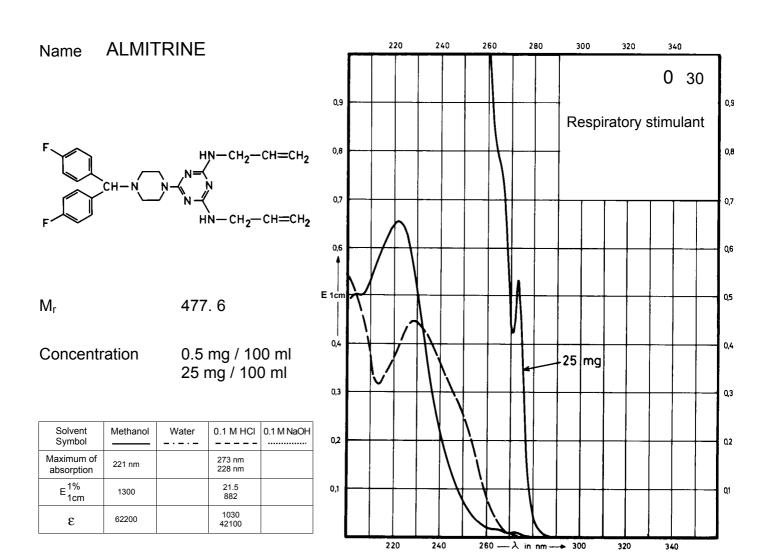

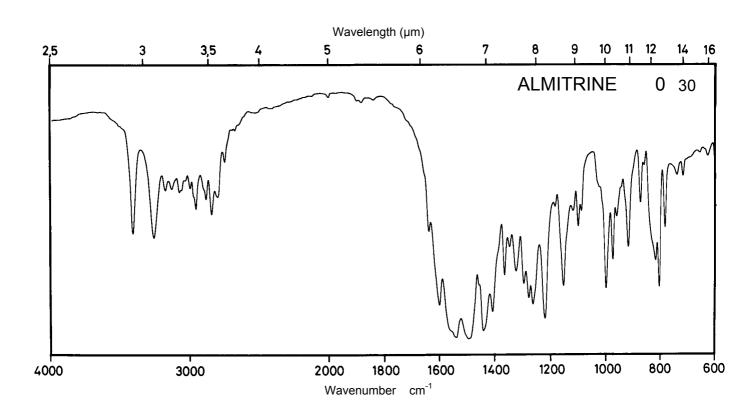



Name AMIPHENAZOLE HYDROCHLORIDE

 M_{r}

227.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 230 nm		261 nm 228 nm	318 nm
E ^{1%} 1cm	417 855		374 935	358
ε	9500 19500		8500 21300	8150

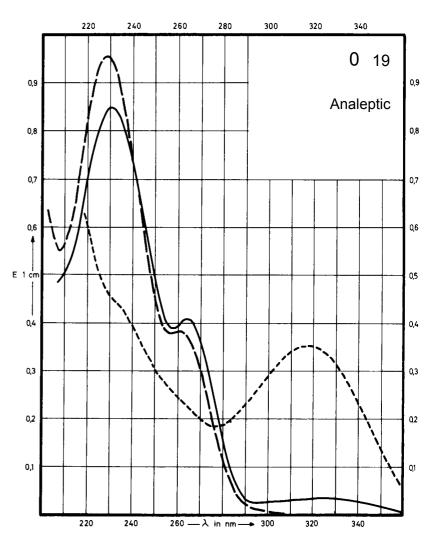

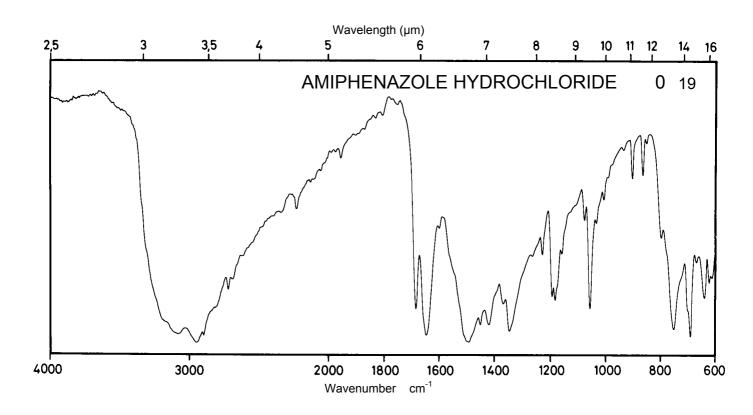



OH NH₂


 M_{r}


109.1

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 232 nm		270 nm 216 nm	291 nm
E ^{1%} 1cm	203 550		180 584	281
ε	2210 6000		1960 6370	3070

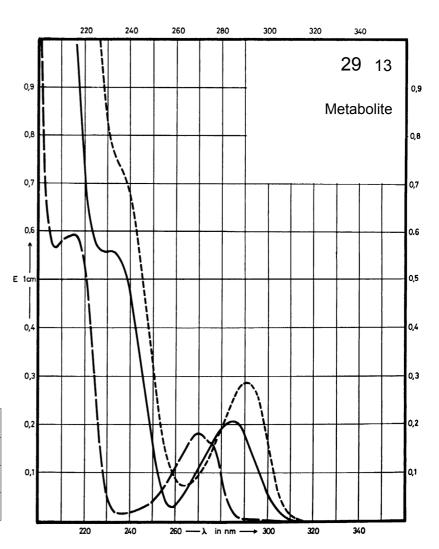

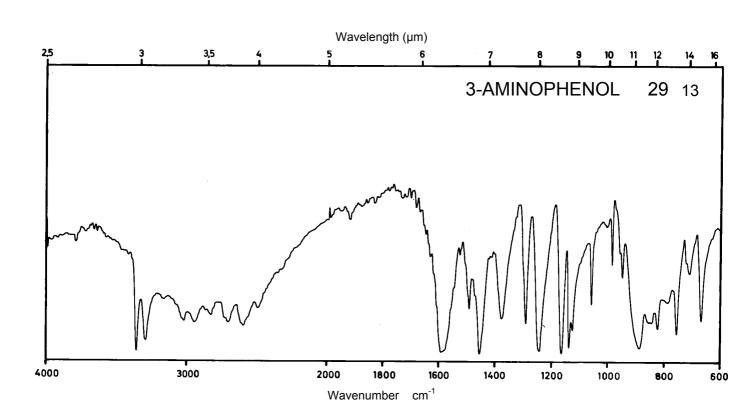





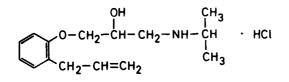


285.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 271 nm		276 nm 270 nm	271 nm
E ^{1%} 1cm	58 64		53 60	62
ε	1650 1840		1520 1710	1770

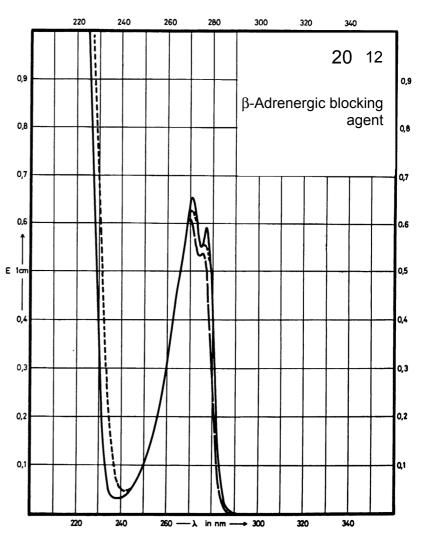

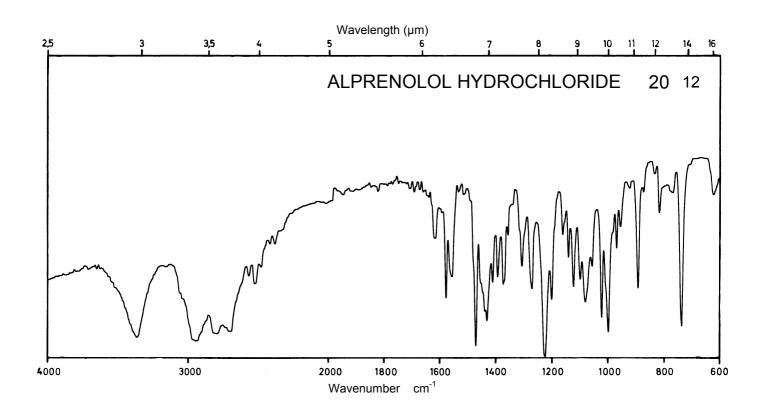





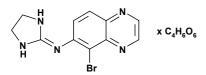



 M_r 442.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	319 nm	319 nm	320 nm	357 nm
	245 nm	247 nm	247 nm	257 nm
E ^{1%}	141	152	154	118
1cm	580	588	569	548
3	6240	6730	6790	5230
	25600	26000	25100	24200

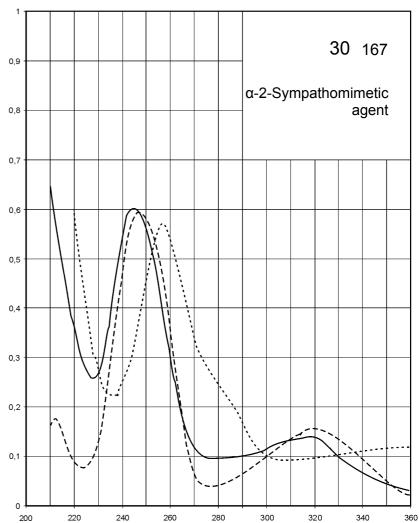


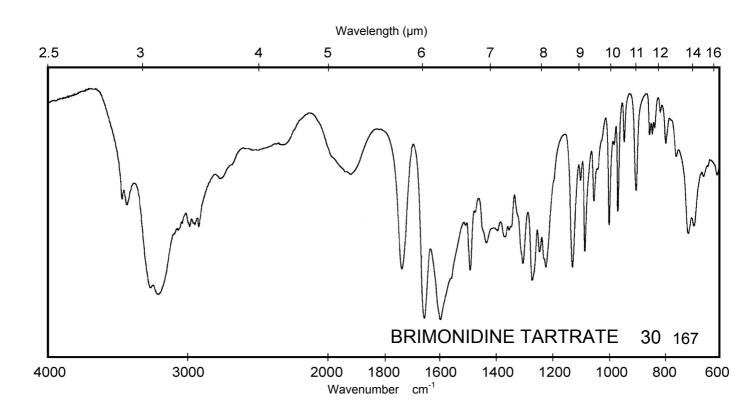

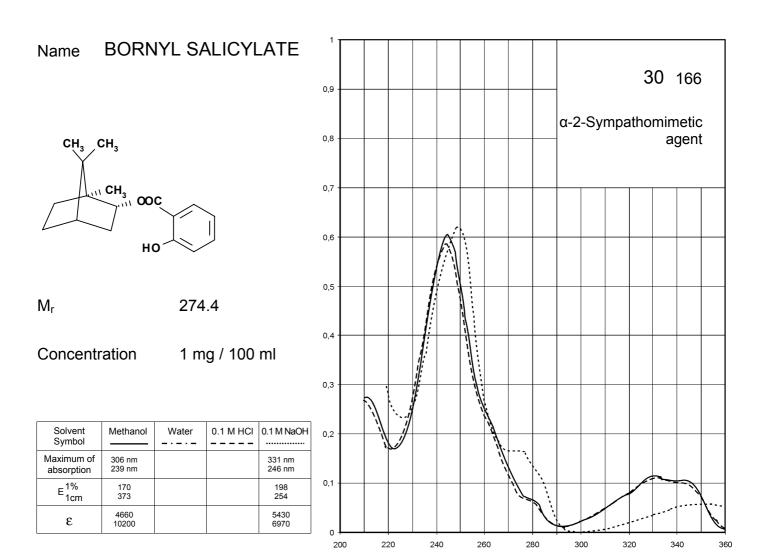

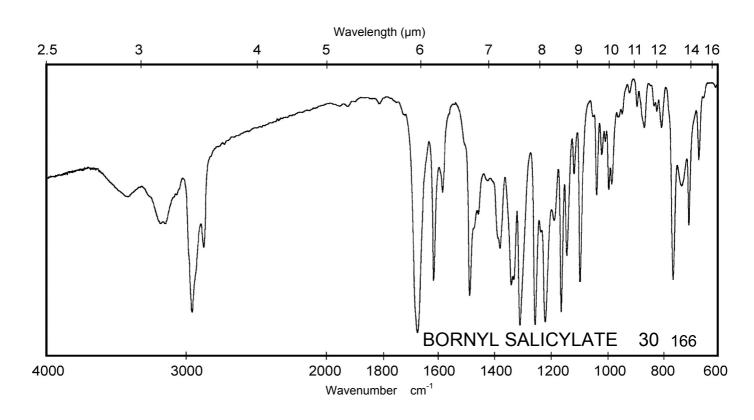



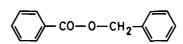



212.2

Concentration


1.26 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 228 nm			
E ^{1%} 1cm	43 683			
ε	910 14500			

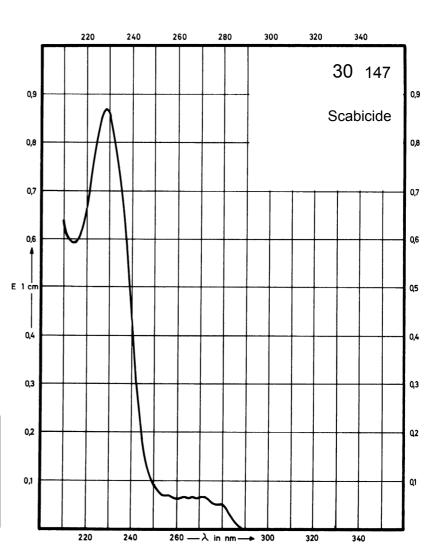

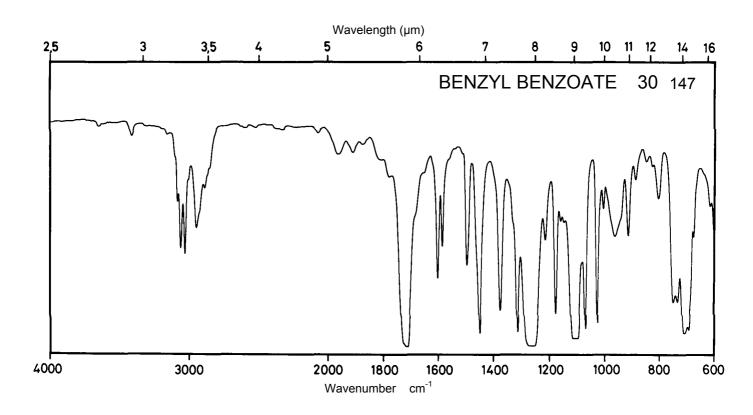



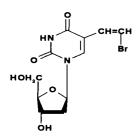



333.1

Concentration

0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	293 nm	291 nm	294 nm	285 nm
	250 nm	250 nm	250 nm	255 nm
E ^{1%}	373	316	335	286
1cm	442	454	456	454
ε	12400	10500	11100	9500
	14700	15100	15200	15100

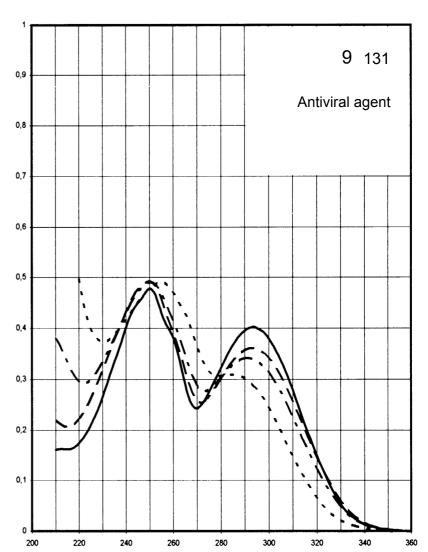

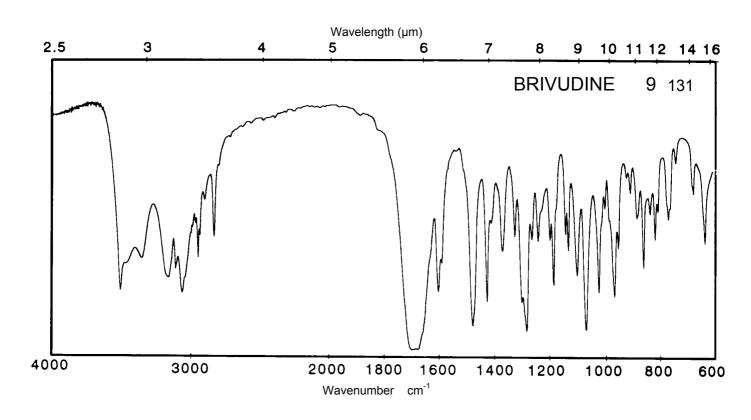



Name BUTYLMETHOXY-PHENOL

OH C(CH₃)₃


 M_r 180.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 227 nm		286 nm 224 nm	305 nm 239 nm
E ^{1%} 1cm	200 326		171 320	215 467
ε	3610 5880		3080 5770	3870 8420

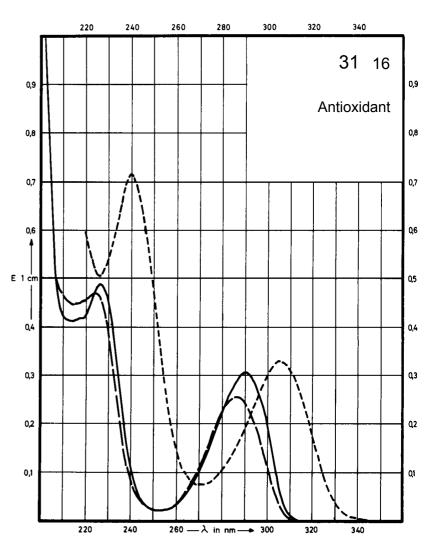

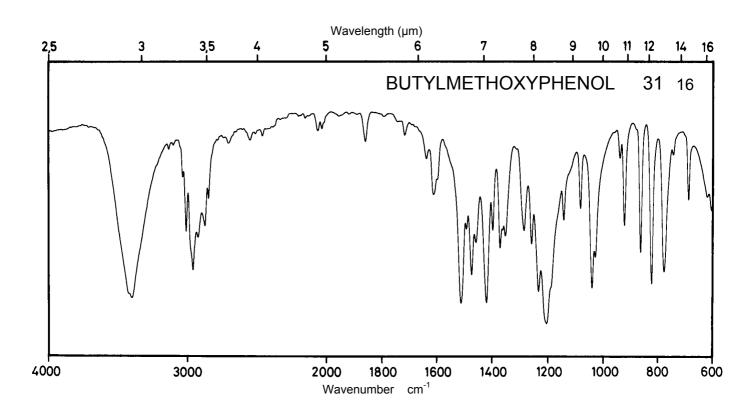



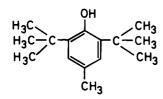



 M_{r}

220.3


Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm			
E ^{1%} _{1cm}	84			
ε	1850			

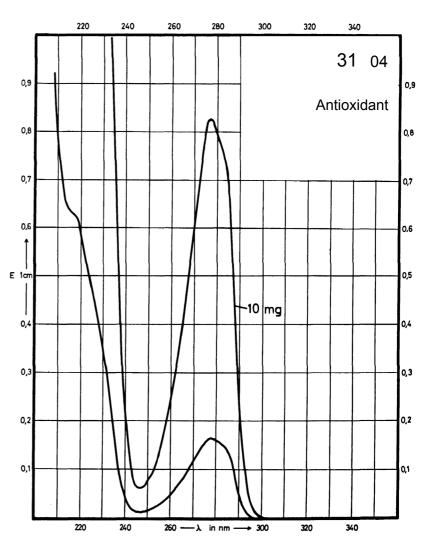

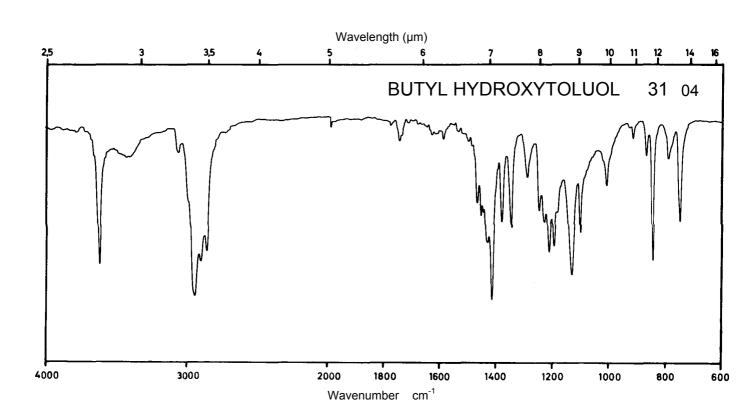





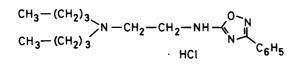


352.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	228 nm		228 nm	
E ^{1%} 1cm	591		552	
3	20850		19500	

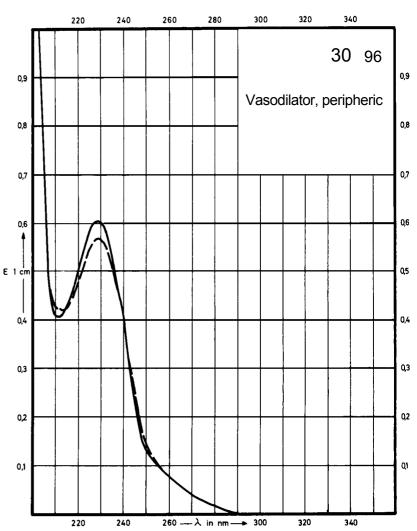

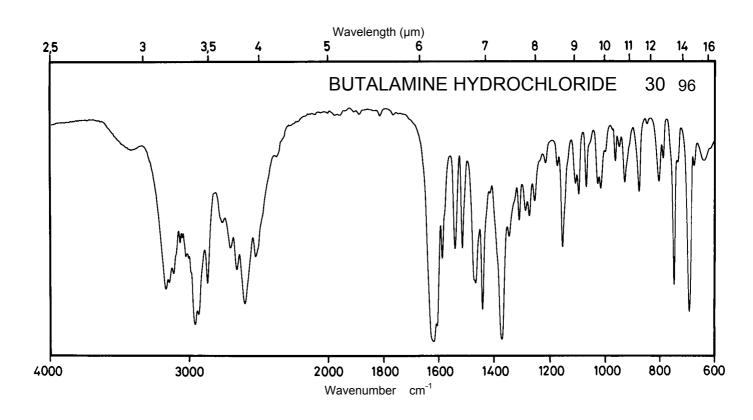





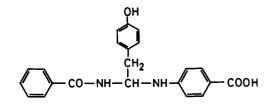

404.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm		270 nm	243 nm
E ^{1%} 1cm	610		532	652
3	24700		21500	26400

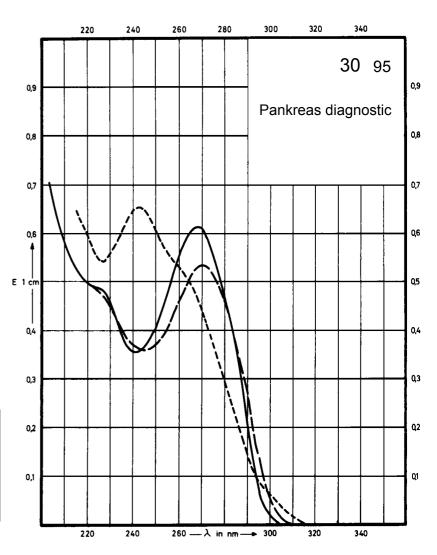

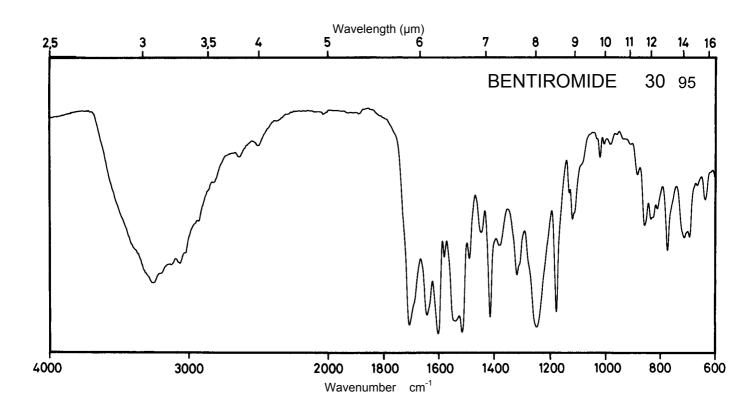



Name BUFLOMEDIL HYDROCHLORIDE

$$\begin{array}{c|c} \mathsf{CO-CH_2-CH_2-CH_2-N} \\ \mathsf{H_3CO} & \mathsf{OCH_3} \\ & \mathsf{HCI} \\ & \mathsf{OCH_3} \end{array}$$


 M_r 343.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm		282 nm	280 nm
E ^{1%} 1cm	146		166	158
3	5000		5700	5400

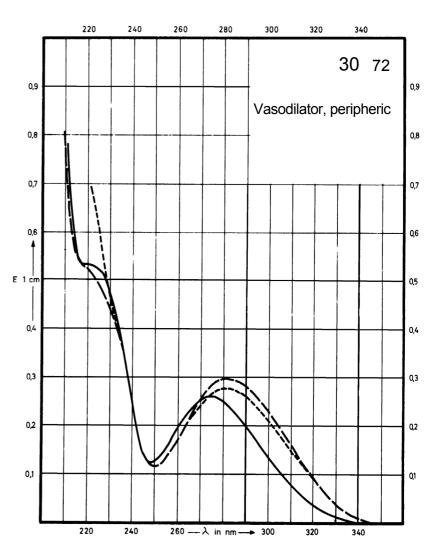

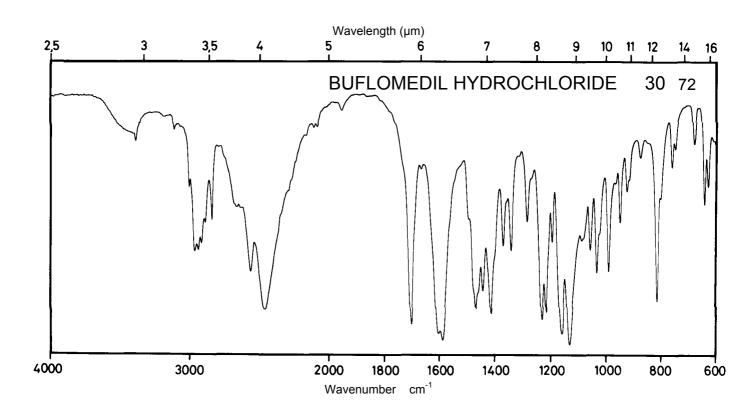





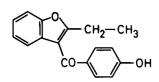


266.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm 231 nm		281 nm 234 nm	353 nm 241 nm
E ^{1%} 1cm	592 680		511 624	931 690
ε	15800 18100		13600 16600	24800 18400

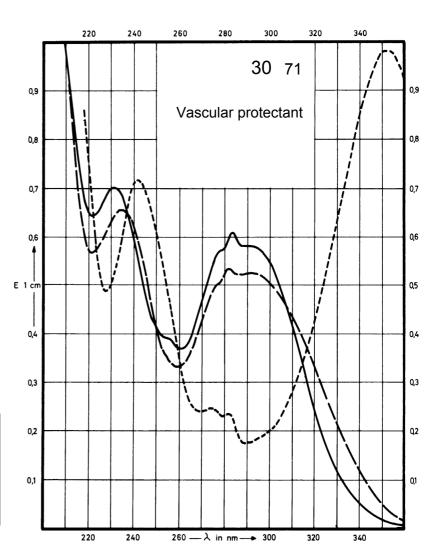

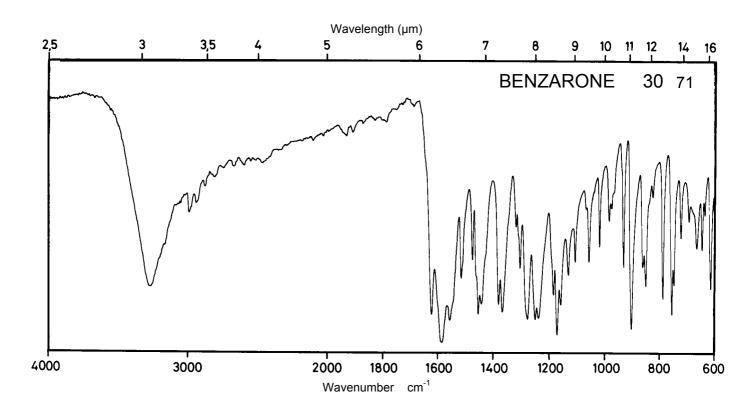

Name BROMOCRIPTINE MESILATE

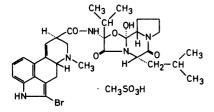


 M_{r}


750.7

Concentration


2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	306 nm	306 nm 240 nm	306 nm	300 nm 239 nm
E ^{1%} 1cm	128	120 273	120	126 325
ε	9640	8990 20470	8990	9500 24360

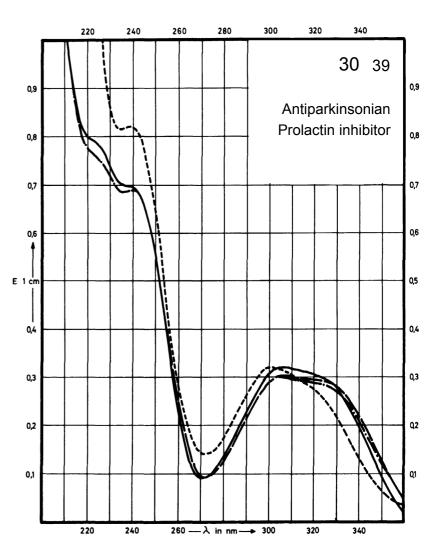

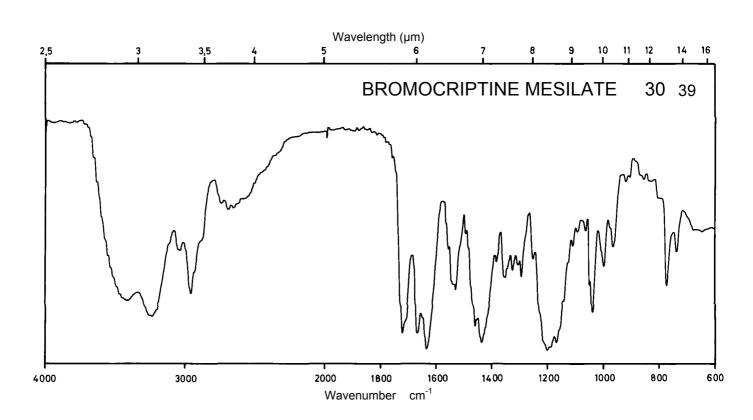





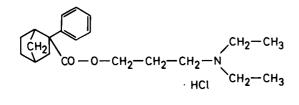


 M_{r} 365.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm 260 nm 254 nm	267 nm 260 nm 254 nm		
E ^{1%} 1cm	4.35 5.70 4.77	4.47 5.95 5.02		
ε	159 209 174	164 218 184		

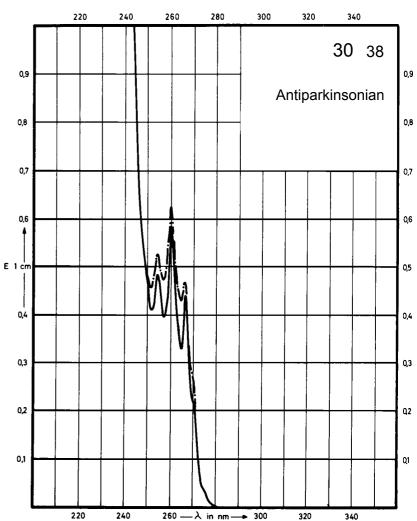

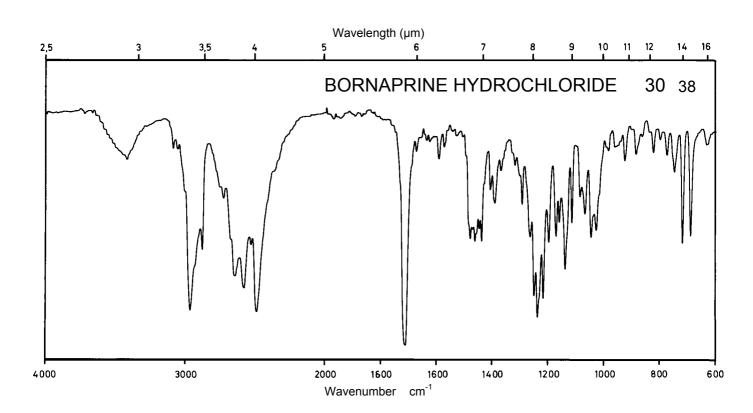




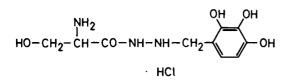


 M_{r} 293.7


Concentration 1 mg / 100 ml 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		272 nm	Decom- position observed
E ^{1%} 1cm	27.4		32.0	
3	805		940	

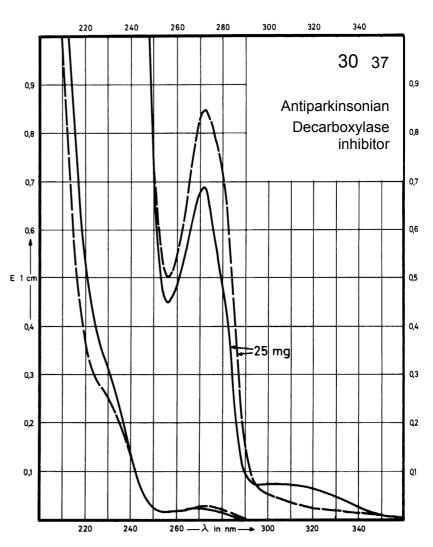

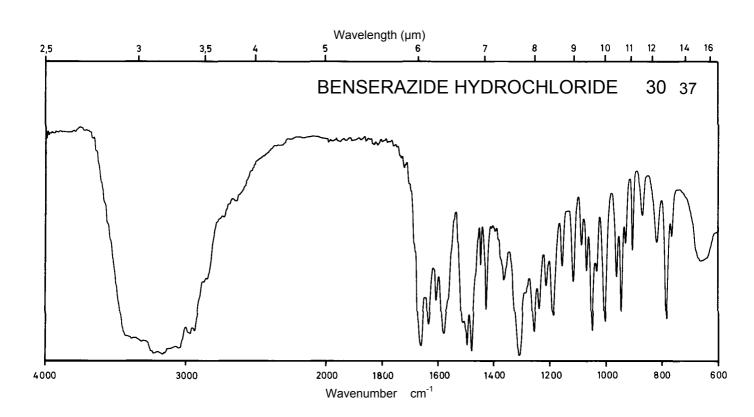



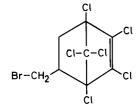



 M_r


393.8

Concentration


3.8 mg / 100 ml

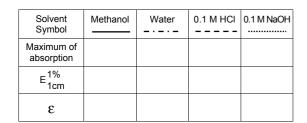

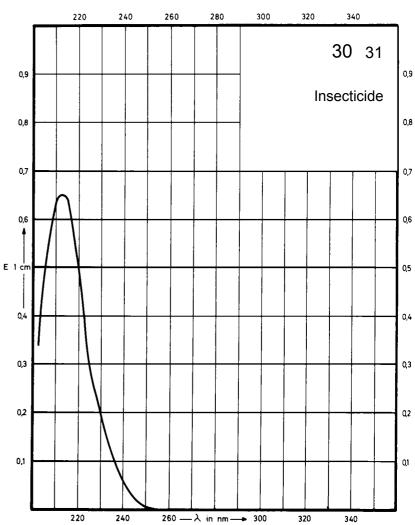

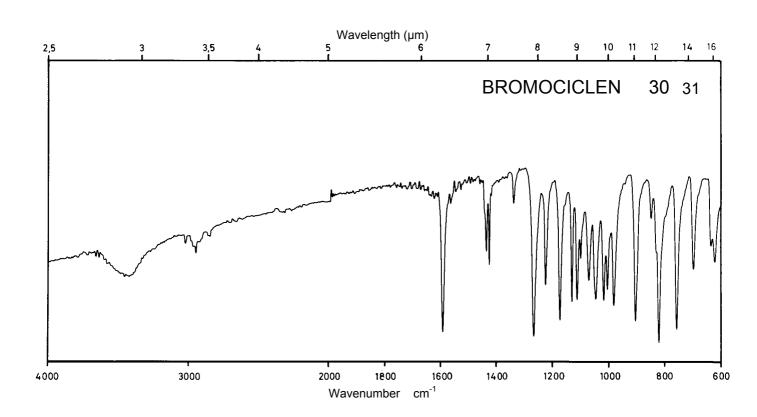



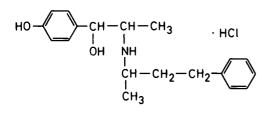



 M_r

335.9

Concentration 2.5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 225 nm		273 nm	291 nm 243 nm
E ^{1%} 1cm	51 300		39	69 390
ε	1710 10070		1320	2330 13080

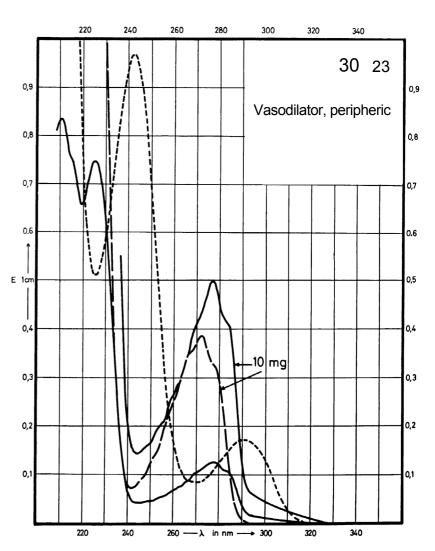

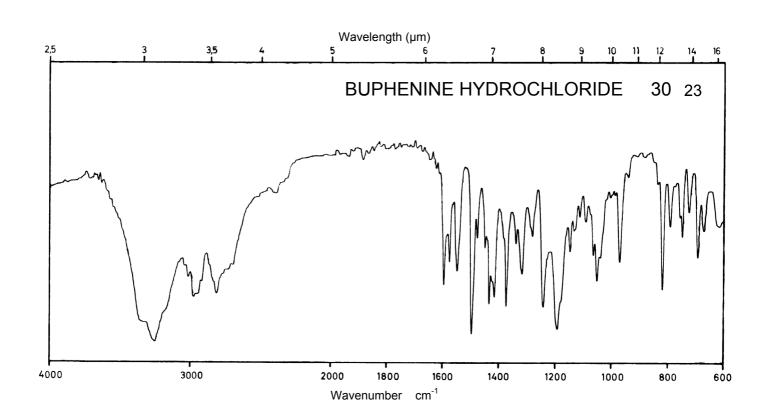


Name BUTYLAMINOBENZOIC ACID


COOH NH C₄H₉

 M_{r} 193.2

Concentration 0.5 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	302 nm 222 nm		226 nm	280 nm
E ^{1%} 1cm	1200 402		673	848
3	23180 7770		13000	16380

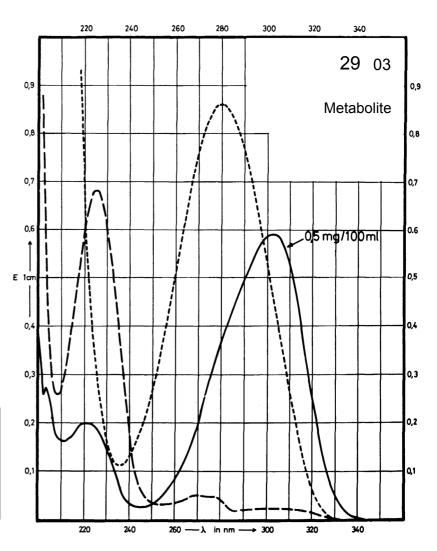

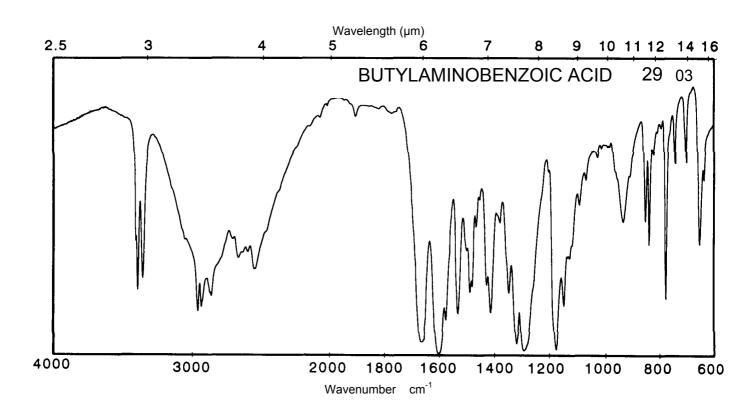




COONa


 M_r

144.1

Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm 222 nm		273 nm 229 nm	
E ^{1%} 1cm	38 624		67 777	
ε	550 9000		970 11200	

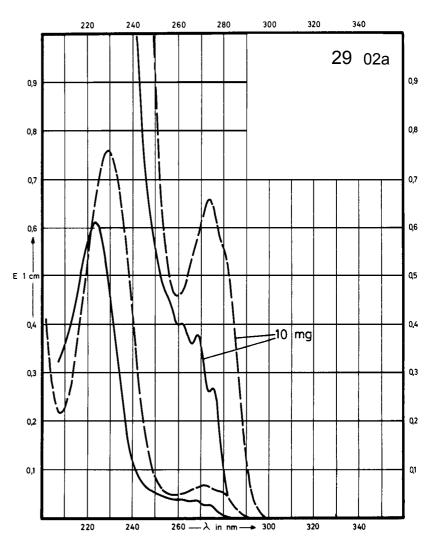

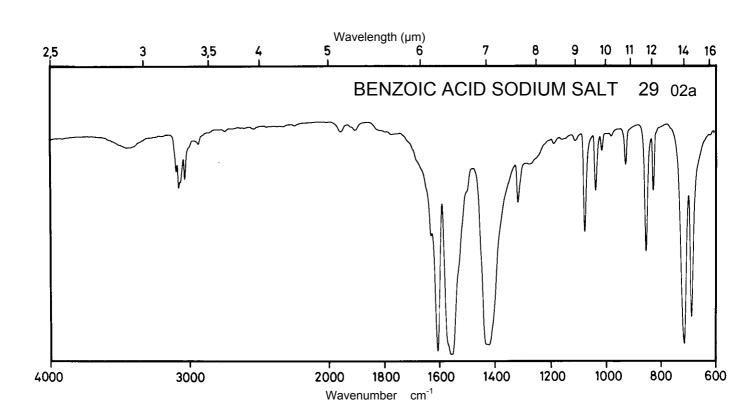





СООН

 M_r 122.1


Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 280 nm 273 nm 227 nm Maximum of 273 nm 230 nm 269 nm absorption 61.4 74.3 E^{1%} 1cm 81.0 929 49.7 931 750 910 11370 ε 610

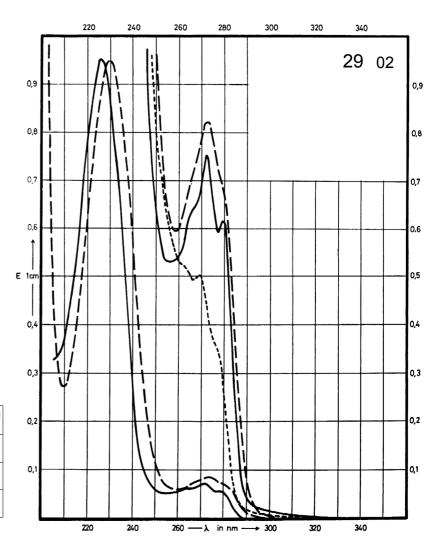

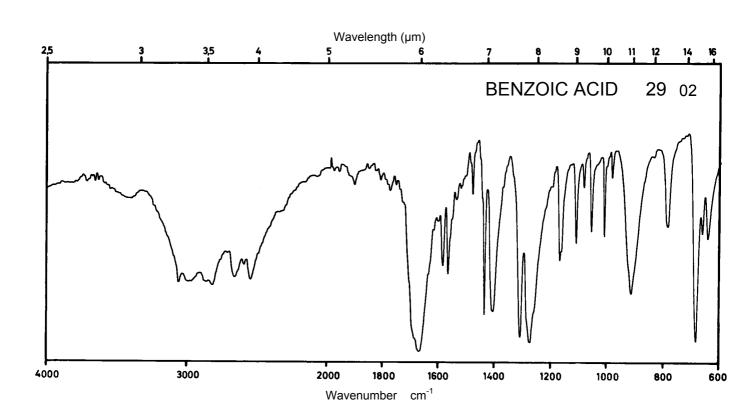





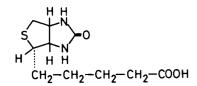


 M_r 244.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

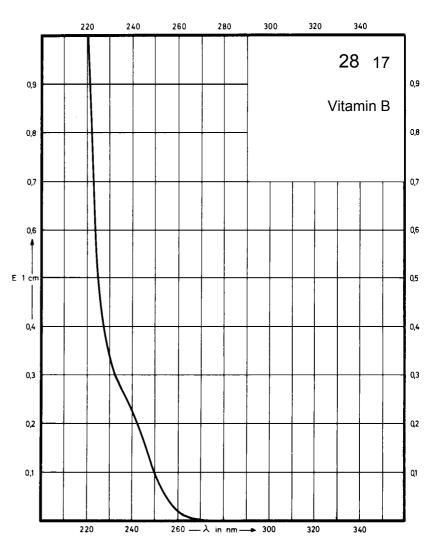

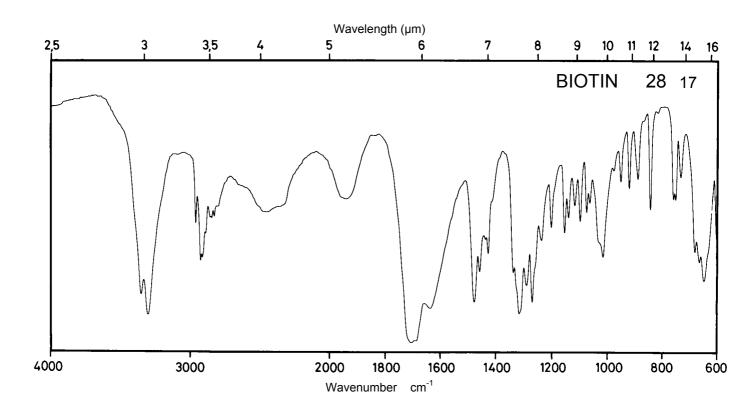





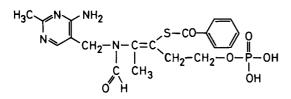


 M_r 466.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm		244 nm	
E ^{1%} 1cm	460		410	
3	21500		19100	

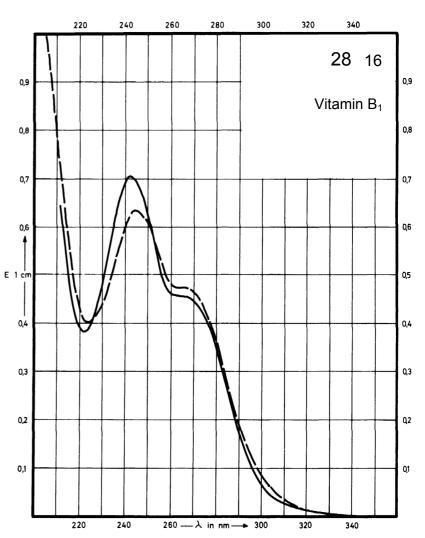

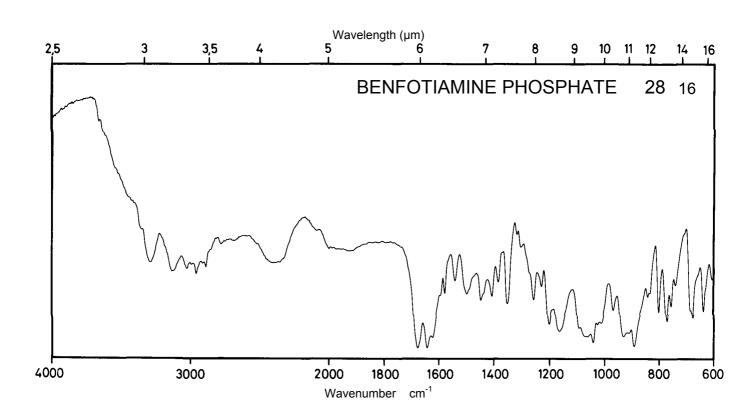




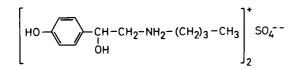


516.6


Concentration 1 mg / 100 ml 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 225 nm		273 nm 222 nm	290 nm 242 nm
E ^{1%} 1cm	62 386		51 348	94 525
3	3200 19940		2630 17980	4860 27120

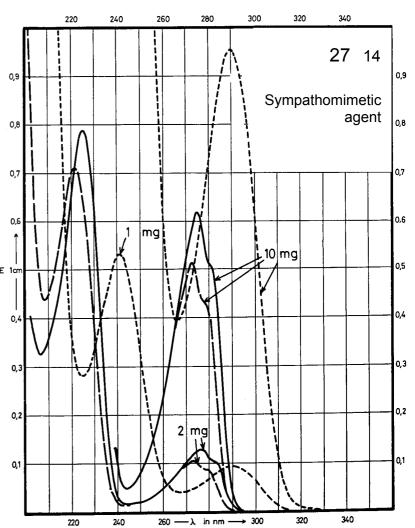

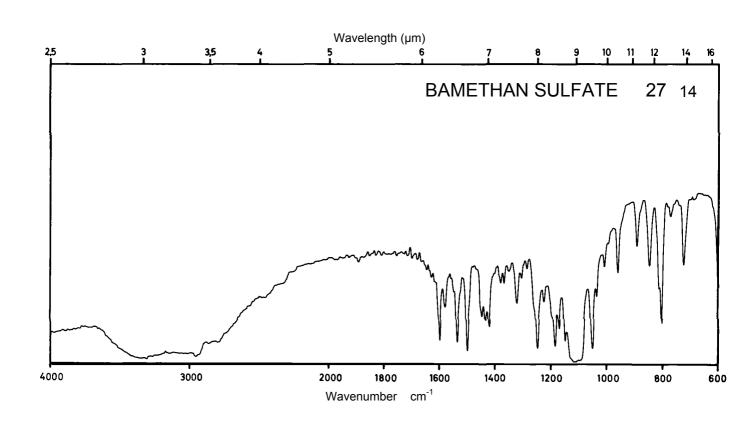


Name BUSPIRONE HYDROCHLORIDE

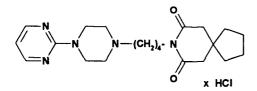


 M_r


422.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	298 nm 236 nm	300 nm 238 nm 213 nm	312 nm 234 nm 213 nm	311 nm 244 nm
E ^{1%} 1cm	56 611	54 530 404	85 529 437	47 527
3	2350 25800	2290 22300 17100	3580 22300 18400	2000 22200

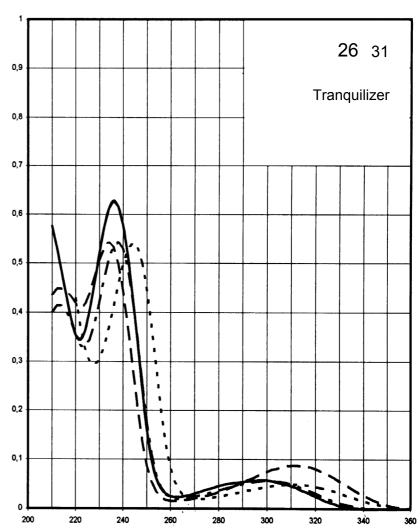

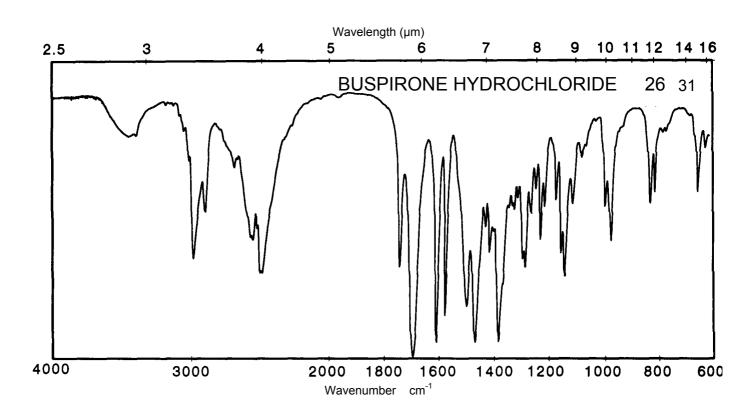



Name BENZOCTAMINE HYDROCHLORIDE

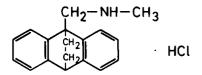

 M_r 285.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 265 nm	272 nm 265 nm	272 nm 265 nm	
E ^{1%} 1cm	46 38	45 38	45 38	
3	1310 1090	1290 1090	1290 1090	

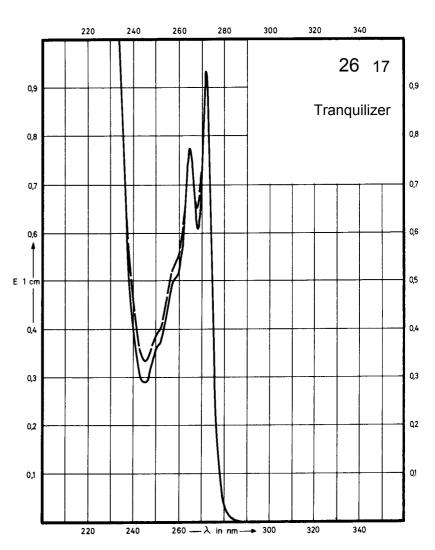

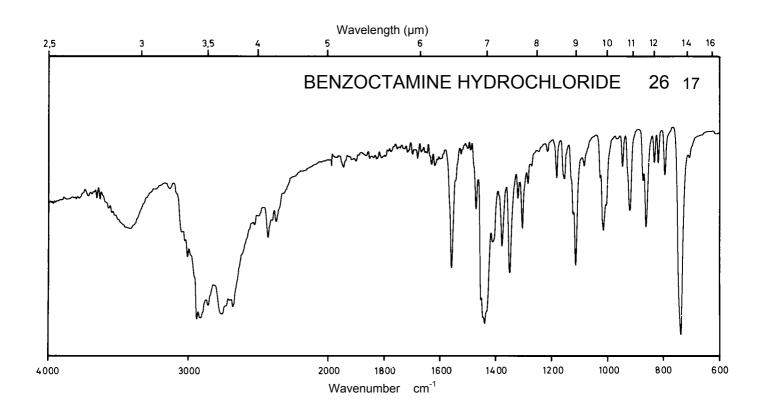


Name BENZBROMARONE


 M_{r}

424.1


Solvent Symbol	Methanol	Water	0.1 M HCI*	0.1 M NaOH
Maximum of absorption	355 nm 237 nm		281 nm 237 nm	355 nm 240 nm
E 1% 1cm	235 562		314 666	513 440
ε	9980 23830		13330 28230	21770 18660

^{* 1}M HCI + Methanol (1 + 9)

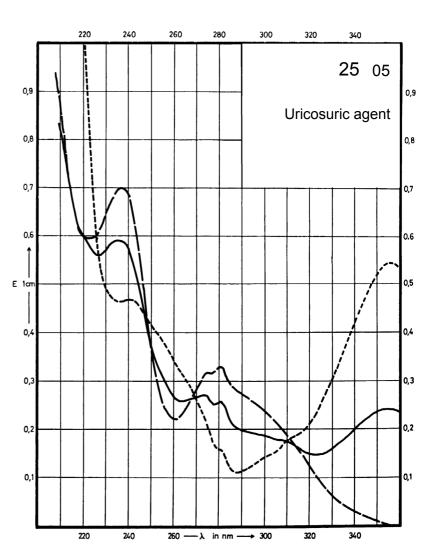

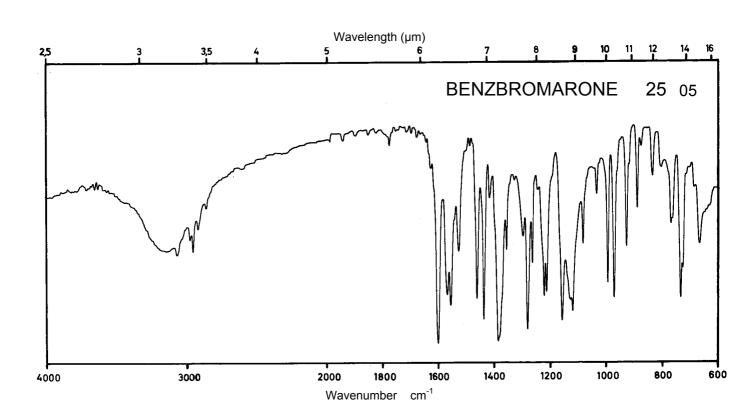




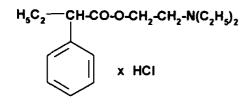


 M_r 299.8


Concentration 1 mg / 100 ml 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 258 nm 252 nm	263 nm 257 nm 252 nm		
E ^{1%} 1cm	4.5 6.2 5.2	4.8 6.5 5.4		
ε	136 186 155	145 194 161		

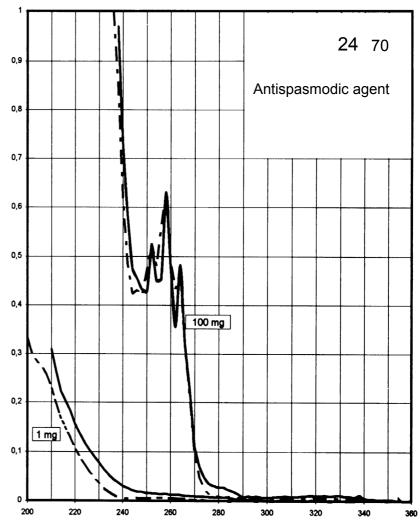

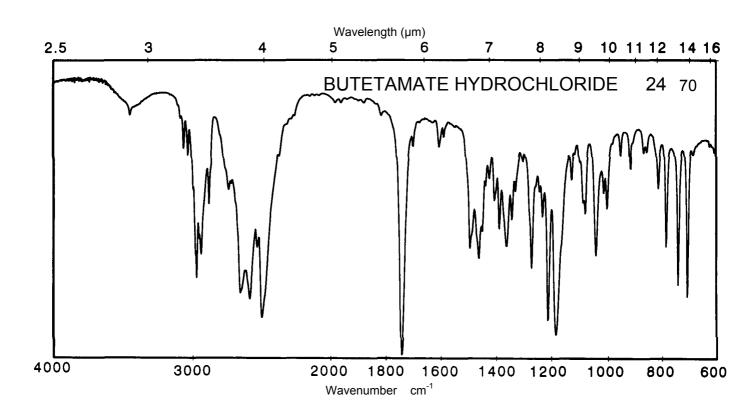





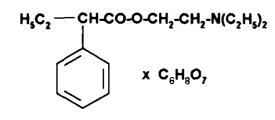


455.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 253 nm	264 nm 258 nm 252 nm		
E ^{1%} 1cm	3.1 4.1 3.4	3.2 4.3 3.7		
ε	142 187	147 197		

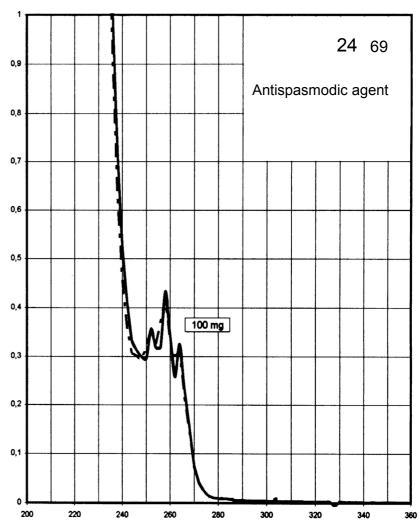

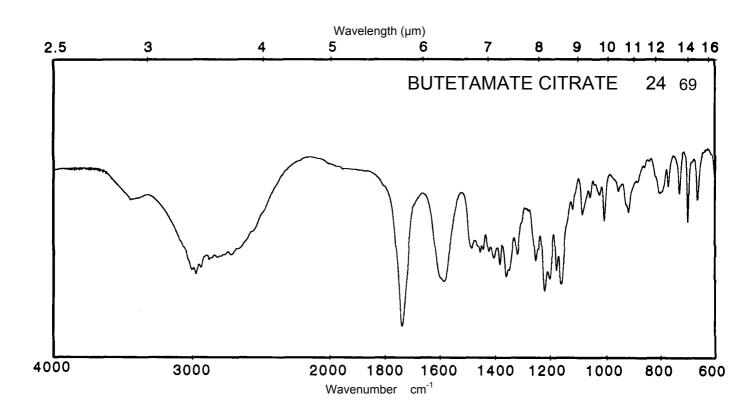


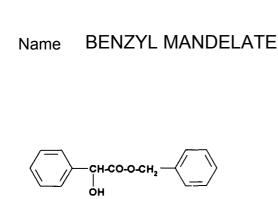


242.3


Concentration 1 mg / 100 ml 12 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm			
E ^{1%} 1cm	13.2 16.5 13.2			
ε	320 400 320			

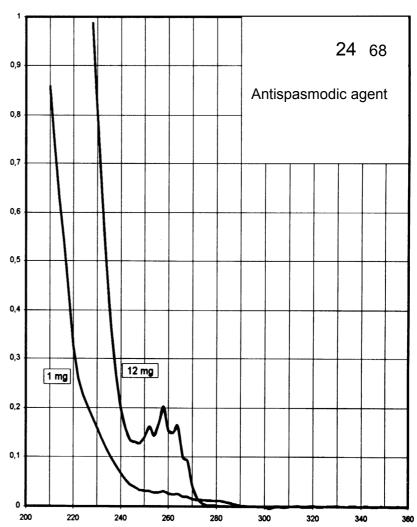

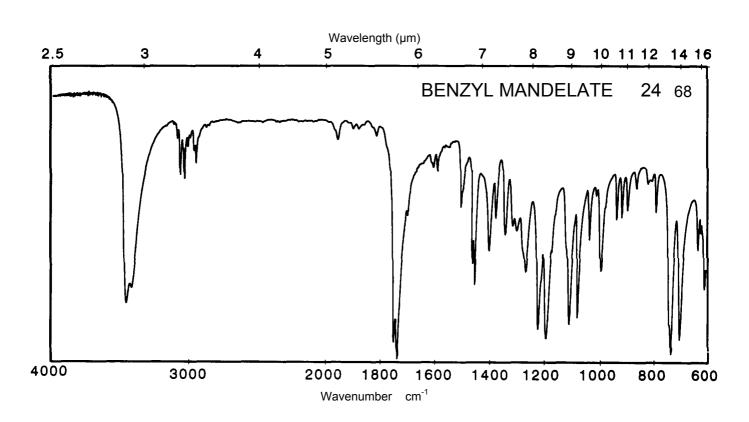




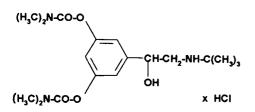


 M_r


403.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm	264 nm	263 nm	
E ^{1%} 1cm	12.3	11.3	11.1	
ε	500	450	450	

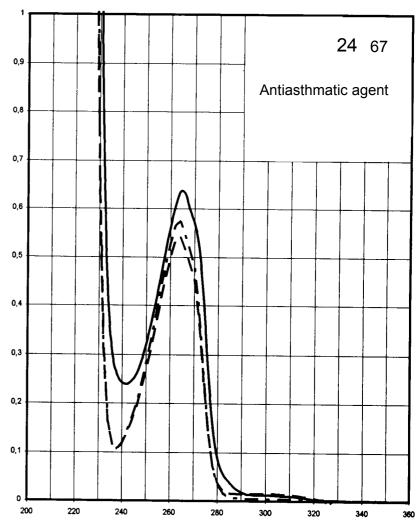

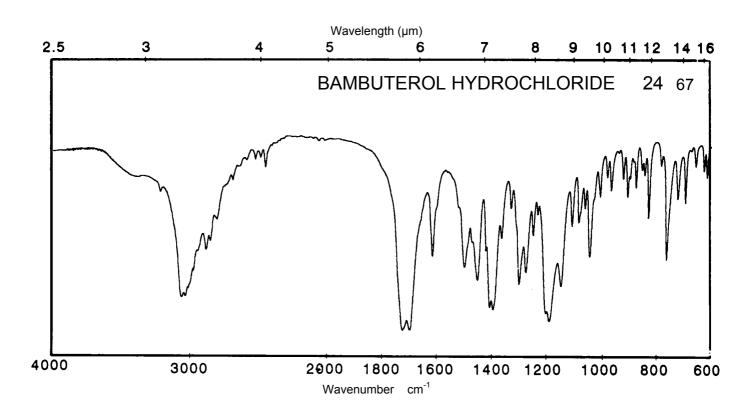





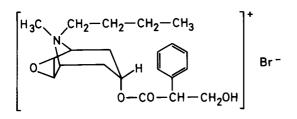


 M_r 440.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm		257 nm	
E ^{1%} 1cm	4.30		4.41	
ε	190		195	

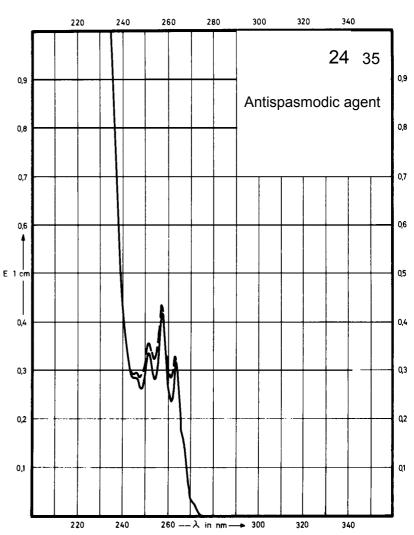

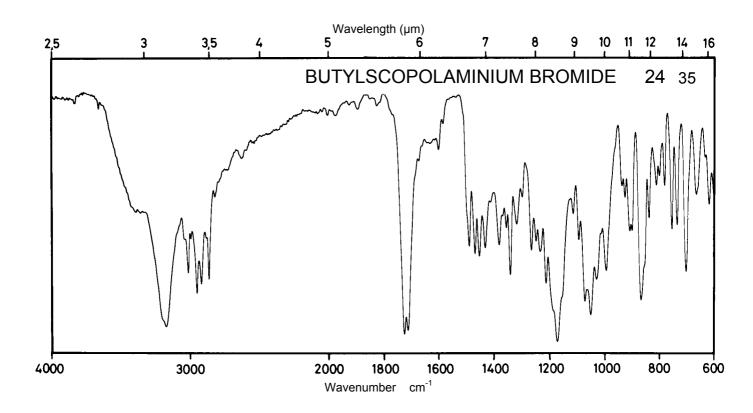




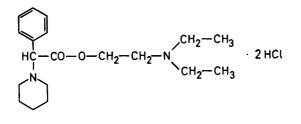


 M_{r} 391.4

Concentration 80 mg / 100 ml 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm 263 nm 258 nm		270 nm 263 nm 258 nm	
E ^{1%} 1cm	13.0 17.1 15.3		13.0 17.4 15.6	
ε	510 670 600		510 680 610	

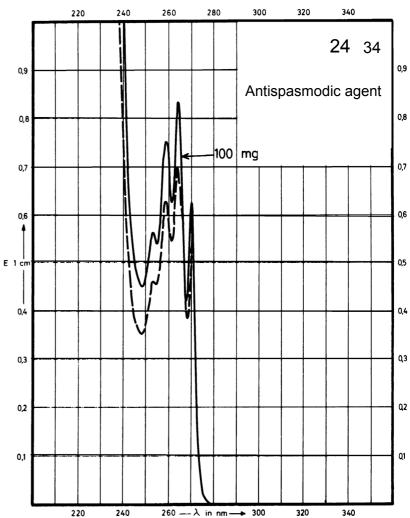

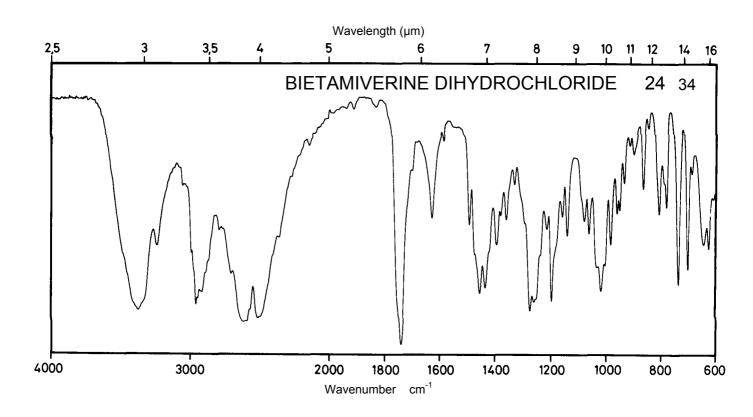


Name BACLOFEN


NH₂-CH₂-CH-CH₂-COOH

 M_{r}


213.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 267 nm 260 nm		274 nm 266 nm 259 nm	275 nm 267 nm 260 nm
E ^{1%} 1cm	9.7 12.4 10.1		8.1 11.4 9.7	11.6 14.8 12.0
3	207 265 215		173 243 207	247 316 257

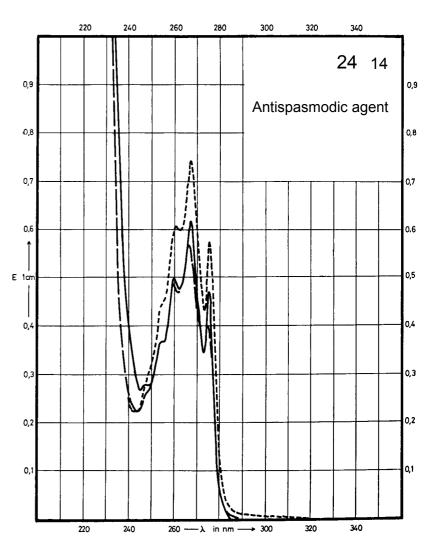

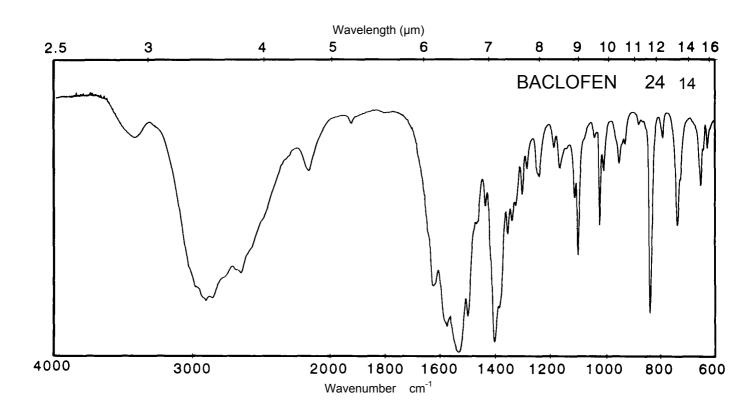




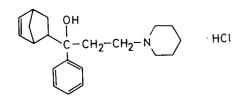


 M_r


347.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm			
E ^{1%} 1cm	4.26 5.56 4.31			
ε	148 194 150			

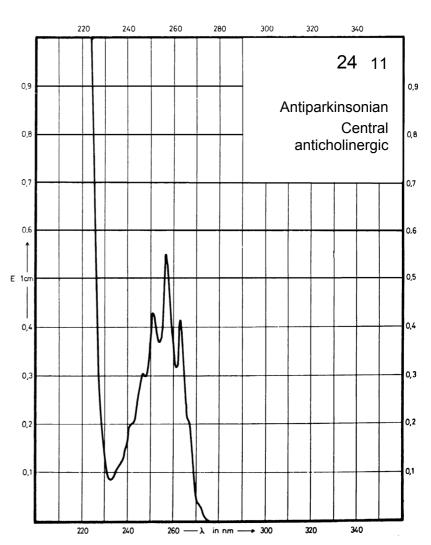

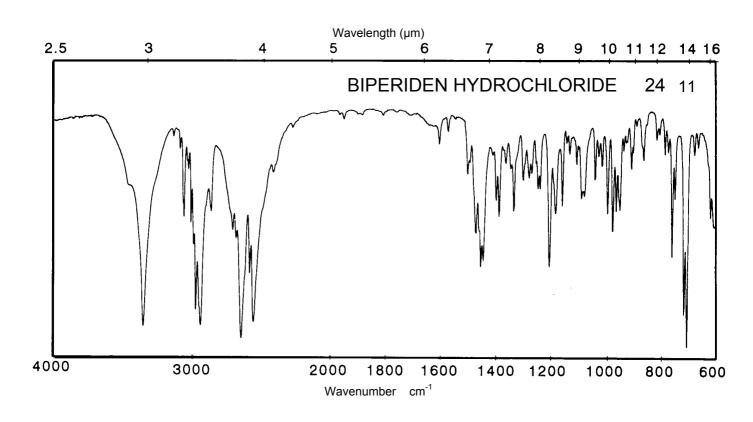



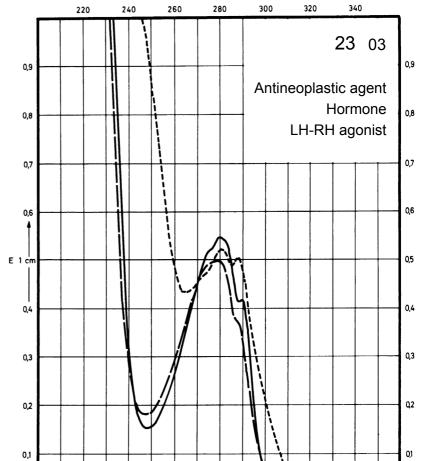


260 — λ in nm → 300

 M_r 1281.4

Concentration 10 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm	279 nm	279 nm	281 nm
E ^{1%} 1cm	54.5	49	49	52
ε	6980	6280	6280	6660

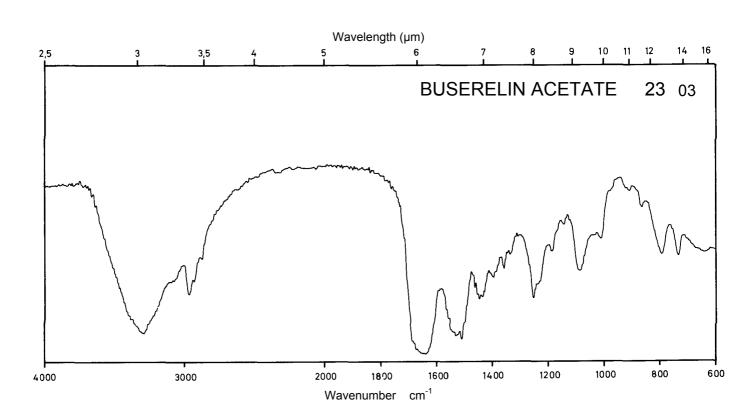

220

240


Name BROMPERIDOL

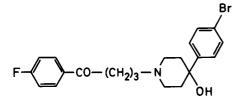

 M_r 420.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm 227 nm 221 nm		247 nm 226 nm	
E ^{1%} 1cm	311 371 373		315 339	
ε	13100 15600 15700		13250 14250	



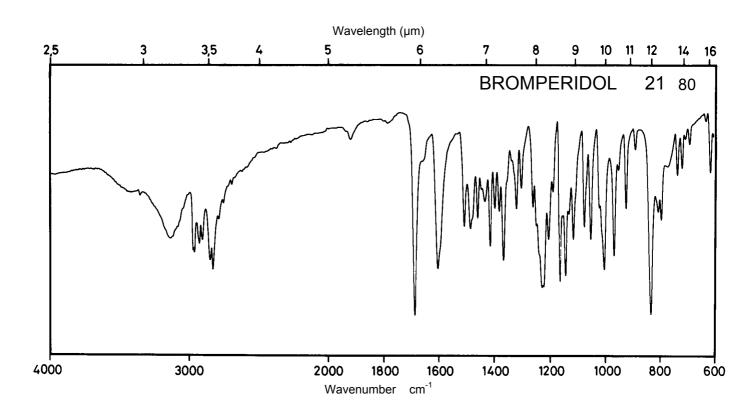



Name BENPERIDOL

 M_{r}

381.5

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	282 nm		279 nm	288 nm
absorption	234 nm		249 nm	248 nm
∟1%	187		204	198
^L 1cm	372		348	445
_	7140		7770	7540
3	14190		13280	16970

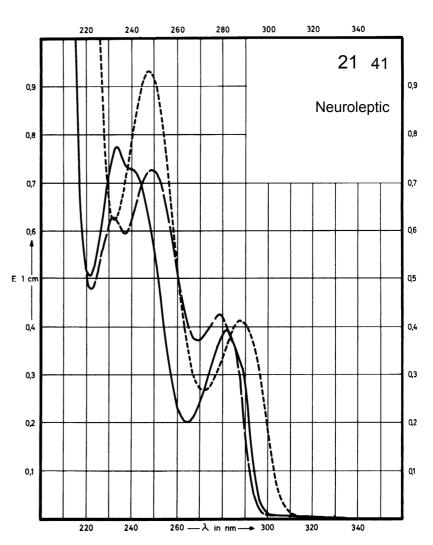

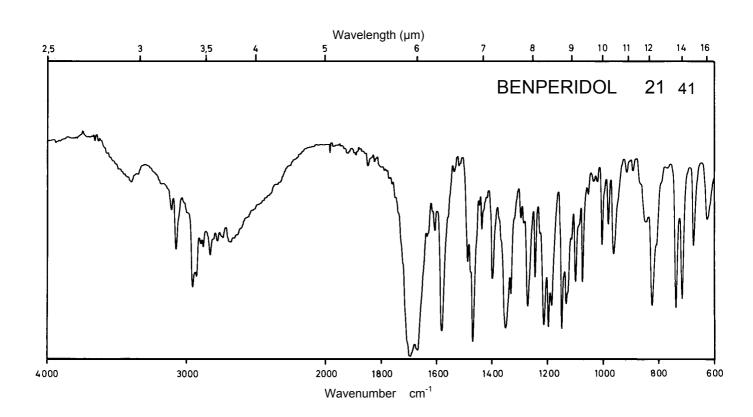



Name BUTAPERAZINE MALEATE

 M_{r}


641.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm 243 nm		277 nm 242 nm	279 nm 242 nm
E ^{1%} 1cm	335 390		288 405	250 353
3	21500 25000		18500 26000	16000 22660

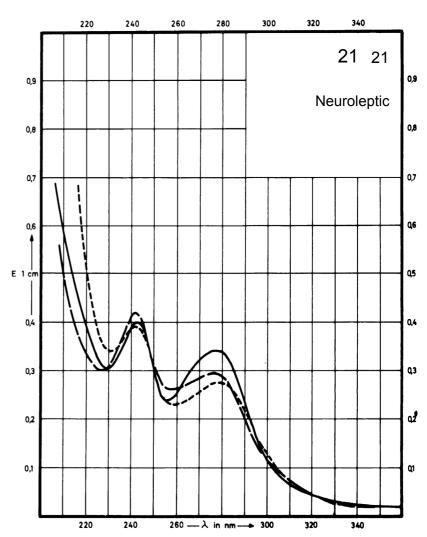

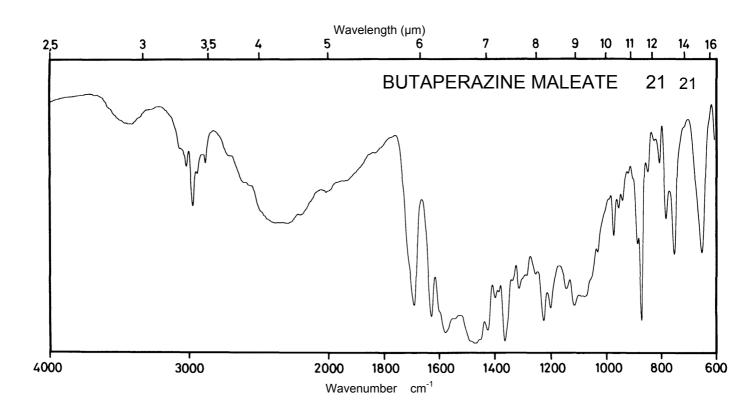




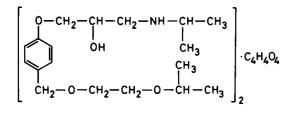


 M_r 767.0


Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm 223 nm		222 nm	271 nm
E ^{1%} 1cm	40 407		409	38
3	3100 31200		31400	2900

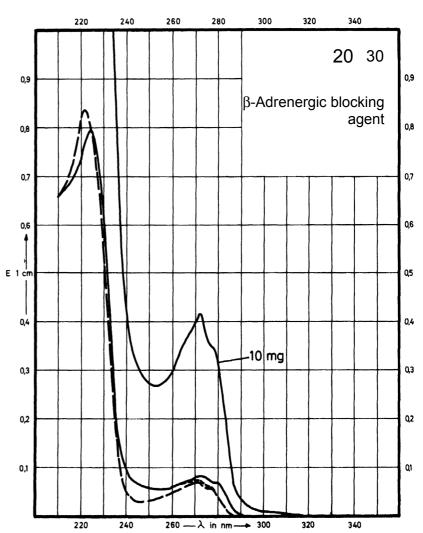

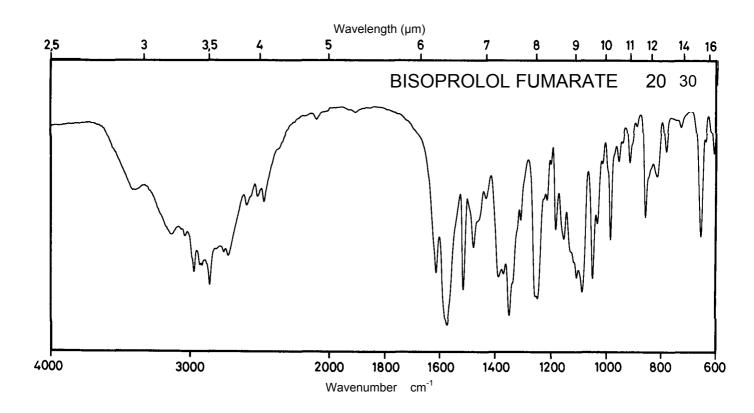


Name BETAXOLOL HYDROCHLORIDE

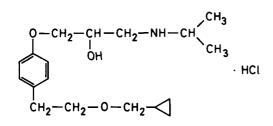


 M_r 343.9


Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 275 nm 222 nm		280 nm 274 nm 221 nm	274 nm
E ^{1%} 1cm	38 46 323		34 41 291	41
ε	1320 1570 11100		1170 1420 10000	1420

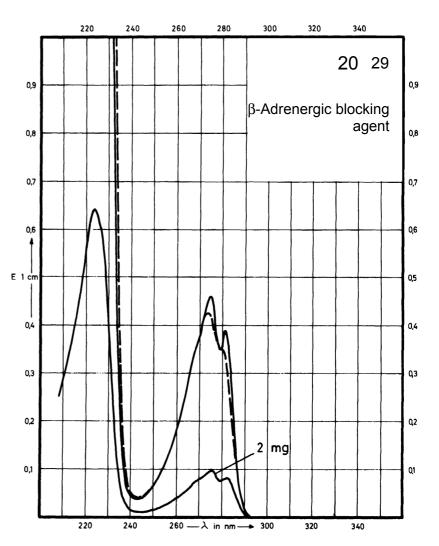

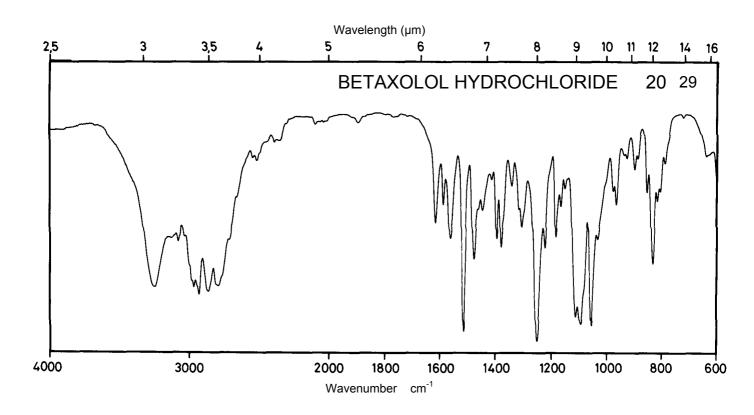





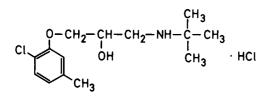


 M_r 308.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm 276 nm		282 nm 275 nm	283 nm 276 nm
E ^{1%} 1cm	55 58		50 54	51 56
ε	1700 1790		1540 1660	1570 1720

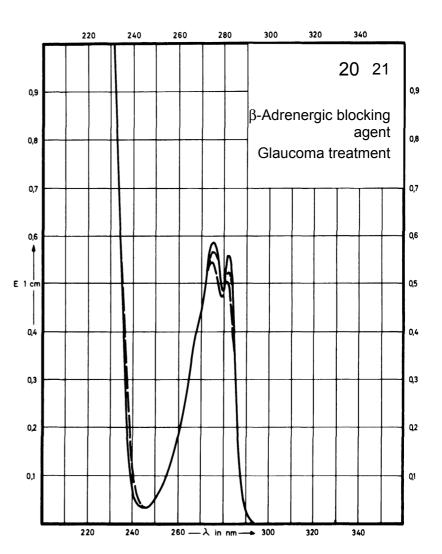

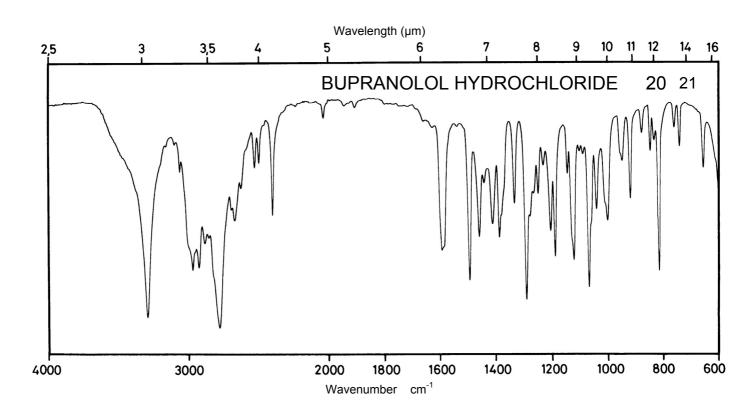


Name BUNITROLOL HYDROCHLORIDE

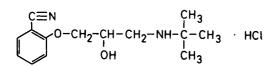


 M_{r}


284.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 231 nm		292 nm 232 nm	294 nm 232 nm
E ^{1%} 1cm	139 322		142 332	149 340
ε	3950 9180		4050 9470	4240 9690

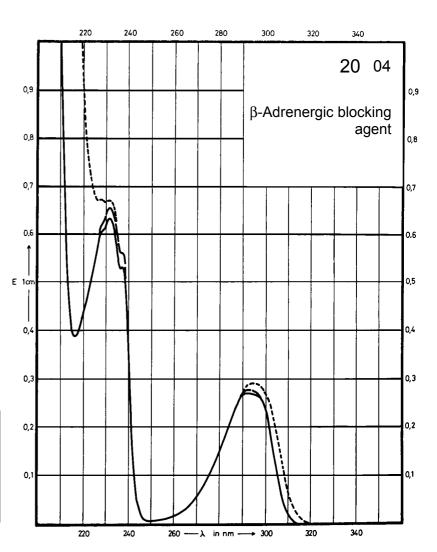

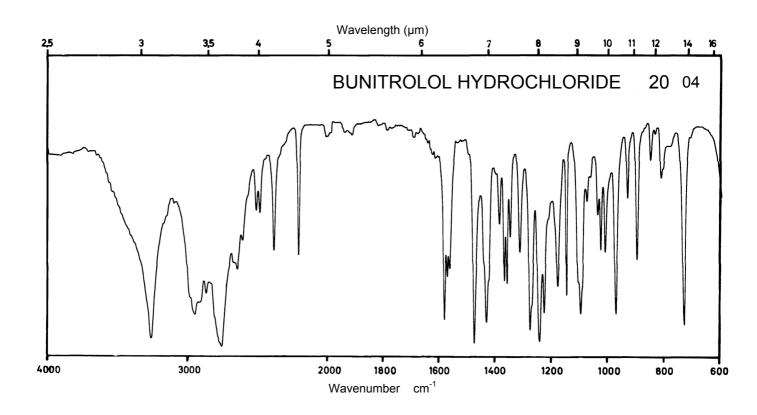





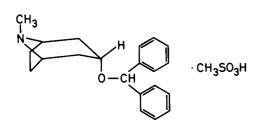

 M_r 403.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm		258 nm	
E ^{1%} 1cm	11.0		11.6	
ε	440		470	

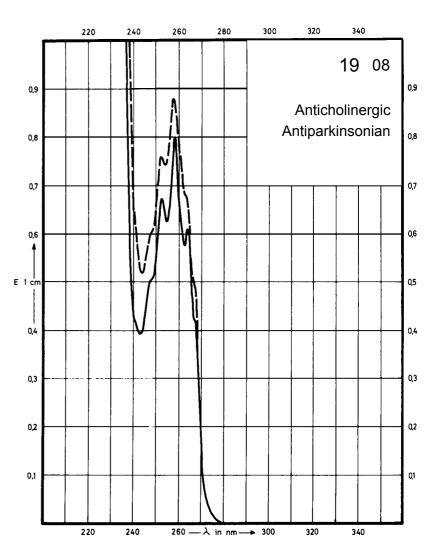

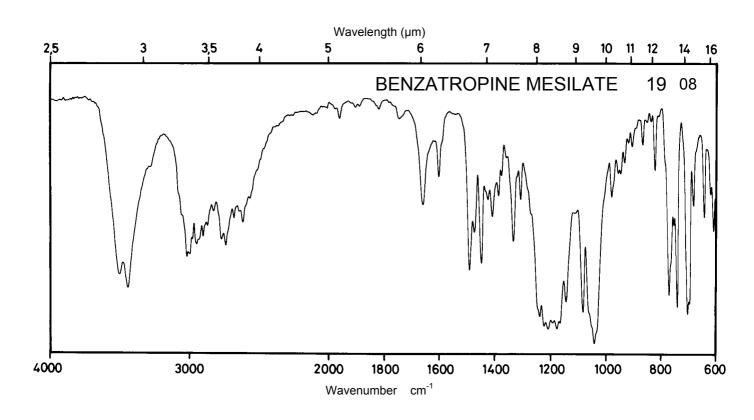





 M_r

361.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	230 nm		226 nm	
E ^{1%} 1cm	578		471	
ε	20920		17030	

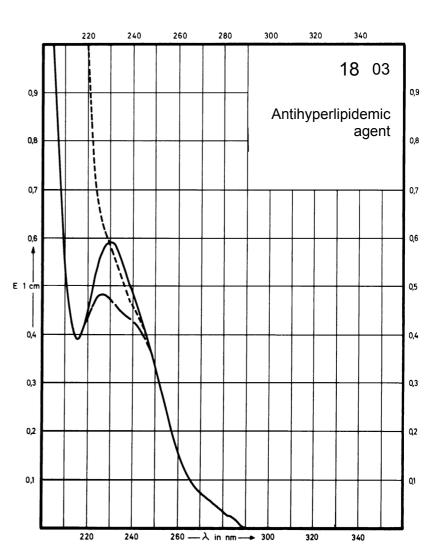

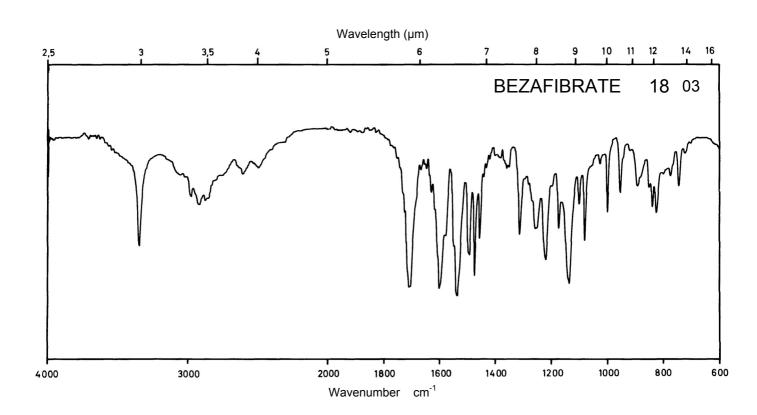



Name BUPIVACAINE HYDROCHLORIDE

 M_{r}


324.9

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of absorption	272 nm 263 nm	271 nm 263 nm	271 nm 263 nm	
E ^{1%} 1cm	8.8 11.5	10.9 13.4	10.9 13.4	
ε	285 373	353 435	353 435	

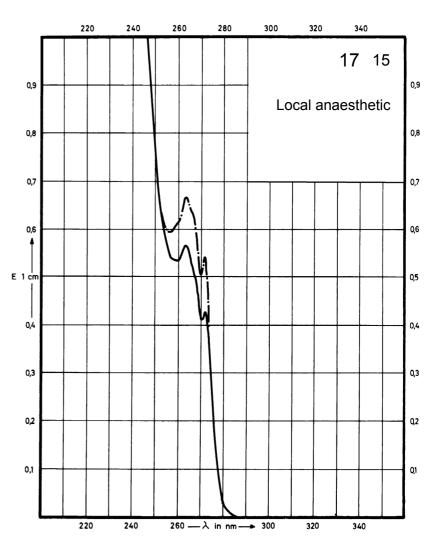

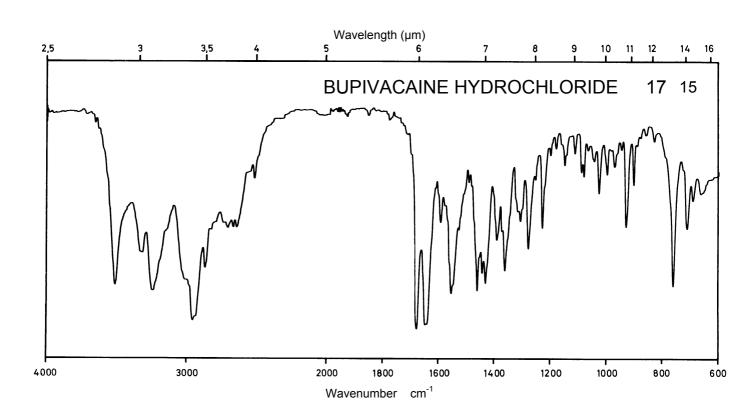





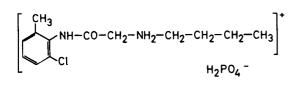


352.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm	267 nm		
E ^{1%} 1cm	9.2	11.2		
3	325	395		

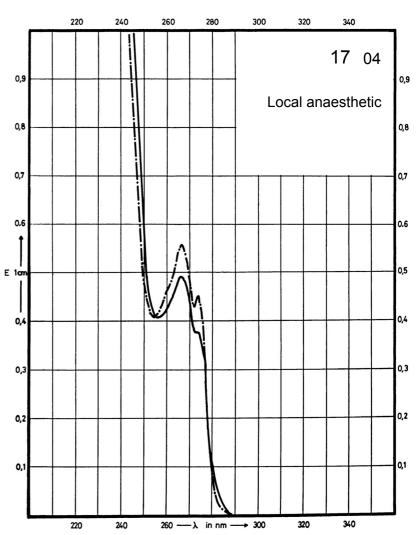

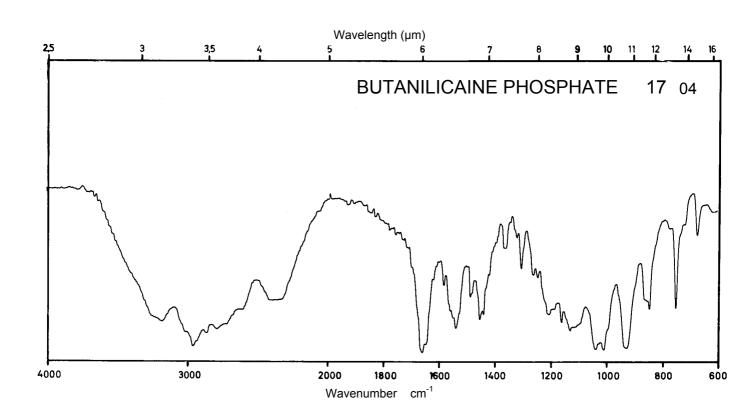





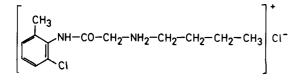


291.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm	267 nm		
E ^{1%} 1cm	11.3	14		
3	330	410		

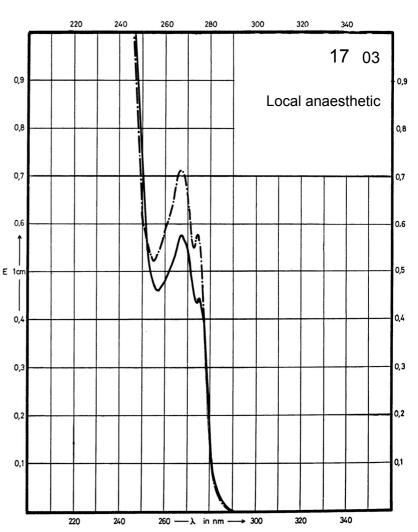

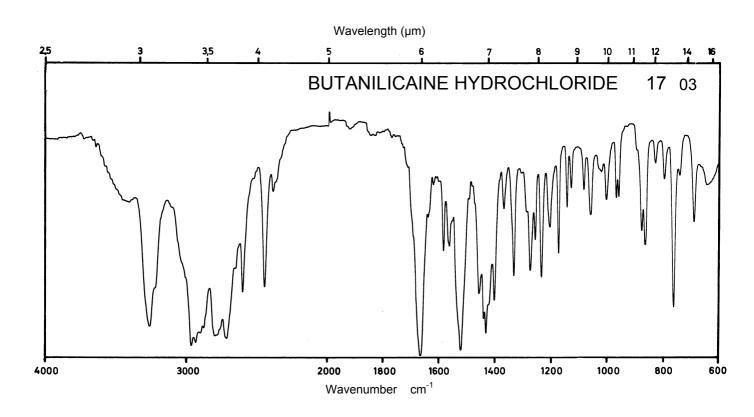





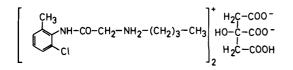

 M_r 701.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm	267 nm		
E ^{1%} 1cm	9.2	11.3		
3	645	790		

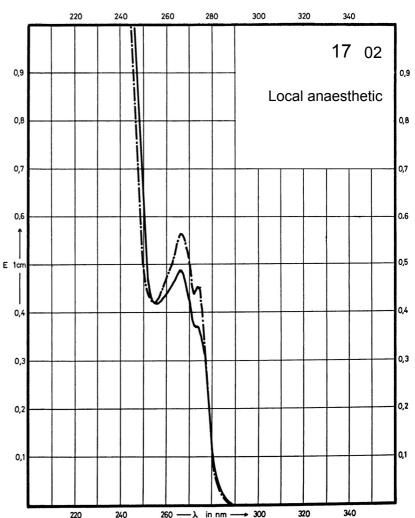

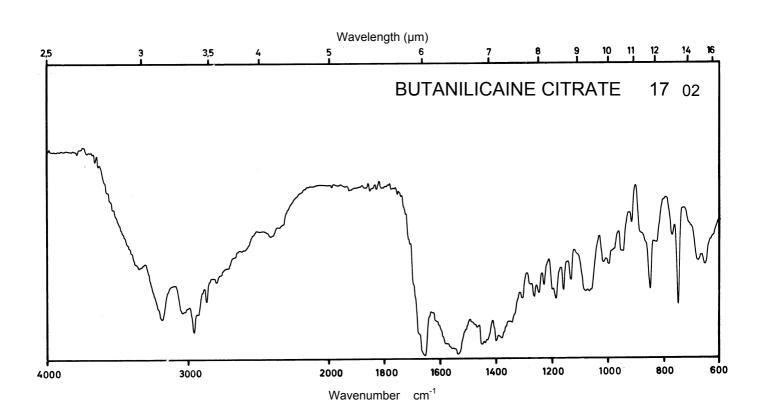







165.2

Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	292 nm 220 nm		270 nm 226 nm	284 nm
E ^{1%} _{1cm}	1246 538		79 770	1002
3	20580 8890		1310 12720	16550

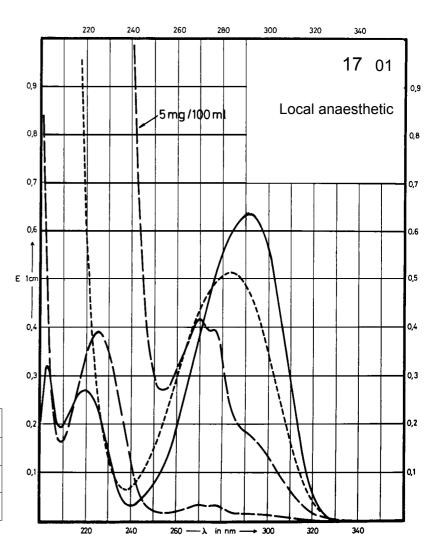

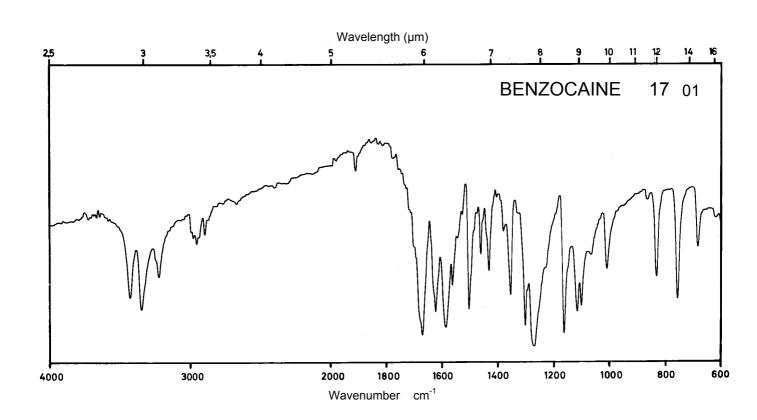



Name BROMOPRIDE

$$\begin{array}{c} \mathsf{CH_2-CH_3} \\ \mathsf{CO-NH-CH_2-CH_2-N-CH_2-CH_3} \\ \mathsf{OCH_3} \\ \mathsf{Br} \\ \mathsf{NH_2} \end{array}$$


 M_r 344.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	309 nm 275 nm		308 nm 273 nm	307 nm 271 nm
E ^{1%} 1cm	381 464		300 334	314 399
3	13100 16000		10300 11500	10800 13700

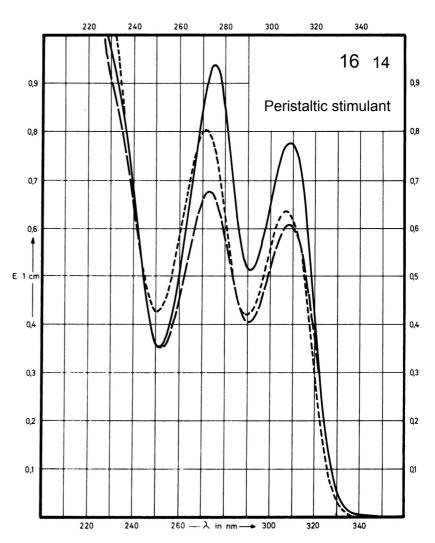

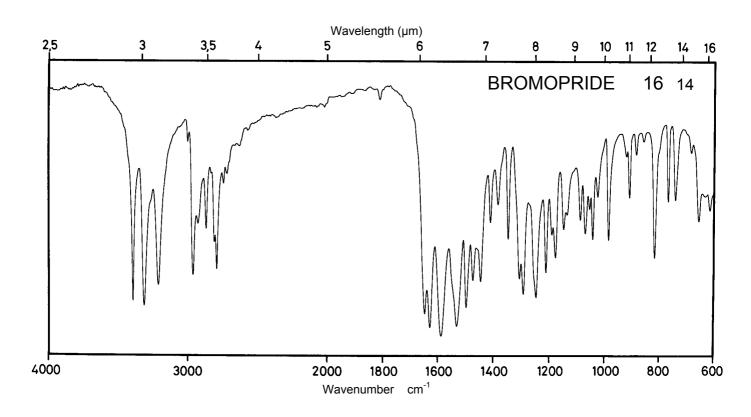




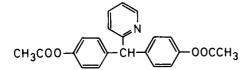


361.8


Concentration 1 mg / 100 ml 2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm		264 nm	246 nm
E ^{1%} 1cm	157		274	618
3	5700		8460	22360

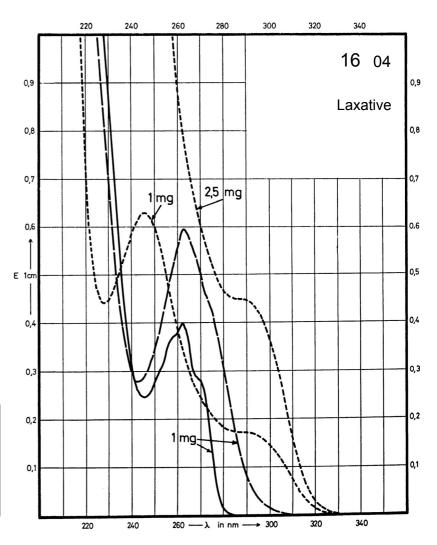

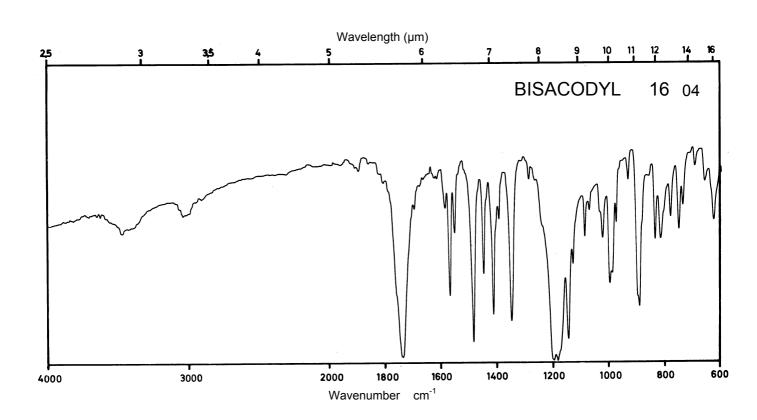



Name BROTIZOLAM

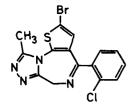

 M_r 393.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm		365 nm 288 nm 255 nm	242 nm
E ^{1%} 1cm	503		73 259 413	518
ε	19800		2870 10200 16260	20400

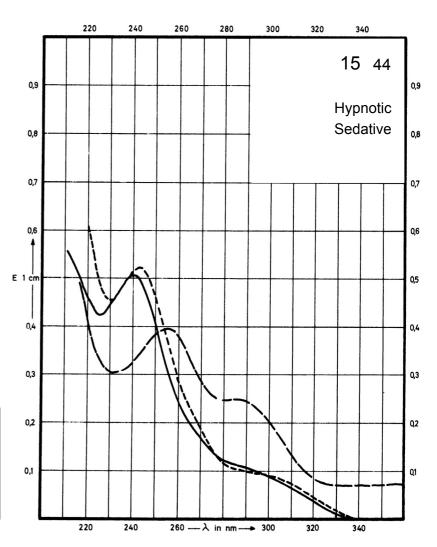

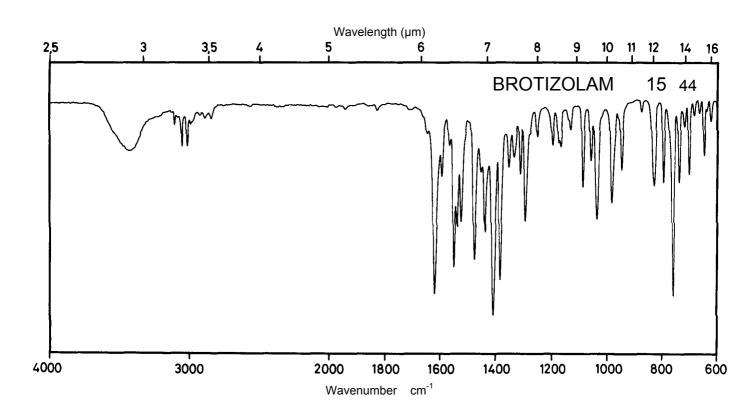

Name BRALLOBARBITAL CALCIUM SALT

$$\begin{bmatrix} O + C & C & C + C + C + C + C + C \\ O + C & C + C + C + C + C + C \\ O & C + C + C + C + C \\ O & C + C + C + C + C \\ O & C + C + C + C \\ O & C + C + C + C \\ O & C + C + C + C \\ O & C + C + C + C \\ O & C + C + C + C \\ O & C$$


 M_{r}


612.3

Concentration

2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				252 nm
E ^{1%} 1cm				189
ε				11600

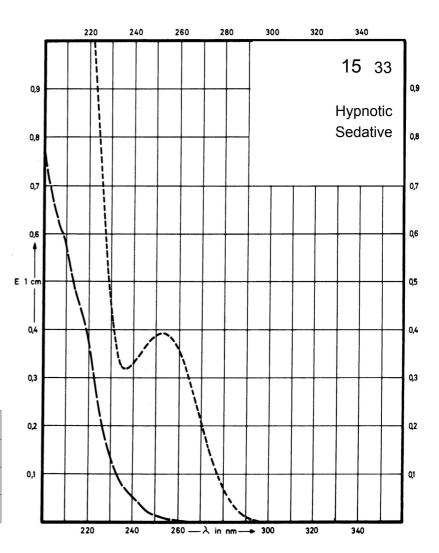

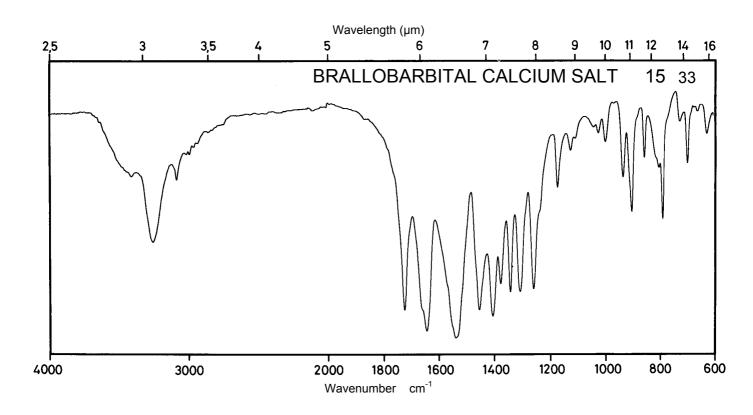





 M_r


223.1

Concentration 5 mg / 100 ml 50 mg / 100 ml

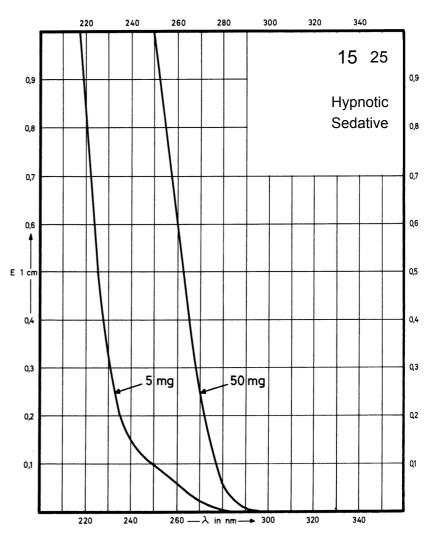

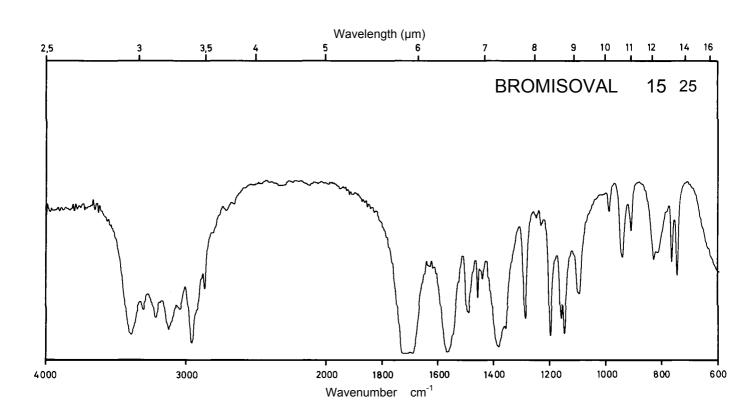





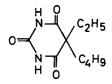


212.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				242 nm
E ^{1%} 1cm				335
3				7110

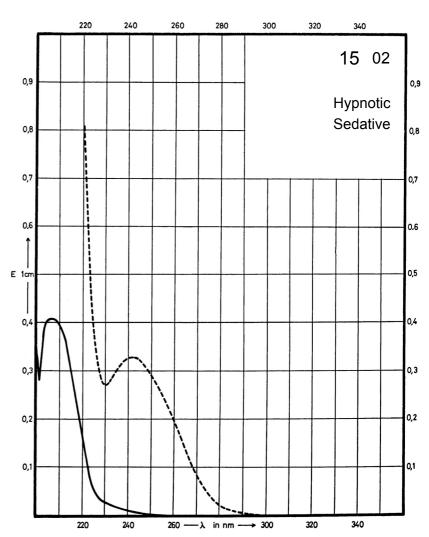

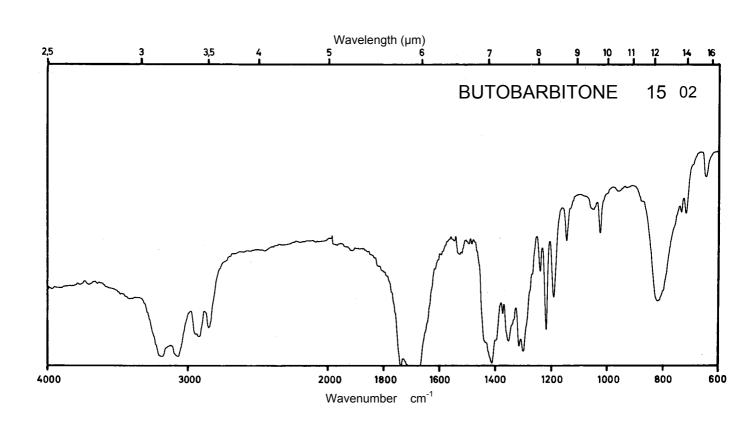





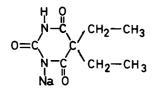


206.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				243 nm
E ^{1%} 1cm				342
3				7050

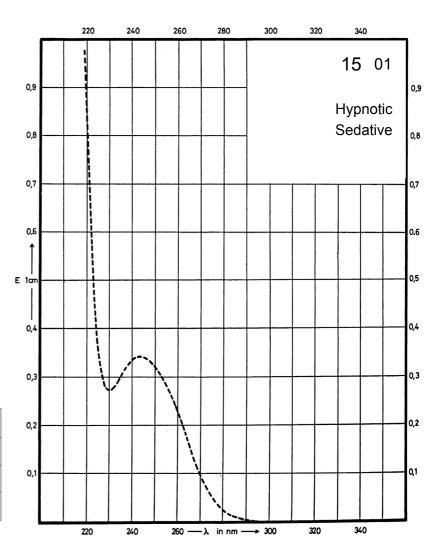

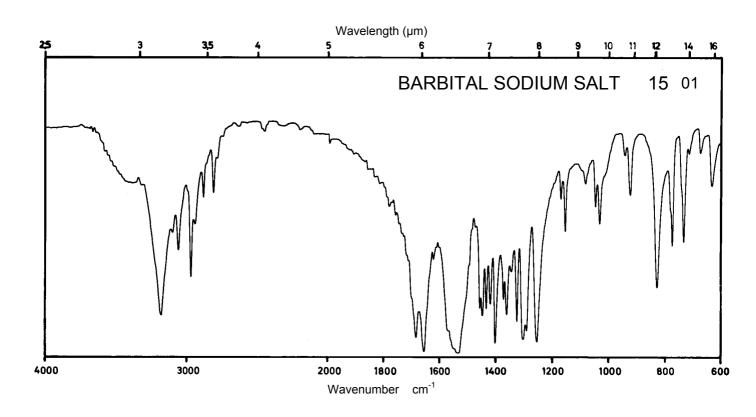





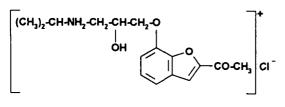


 M_r 327.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	290 nm	294 nm	294 nm	295 nm
	240 nm	240 nm	240 nm	241 nm
E ^{1%}	549	555	555	544
1cm	413	424	423	428
3	18000	18200	18200	17800
	13500	13900	13900	14000

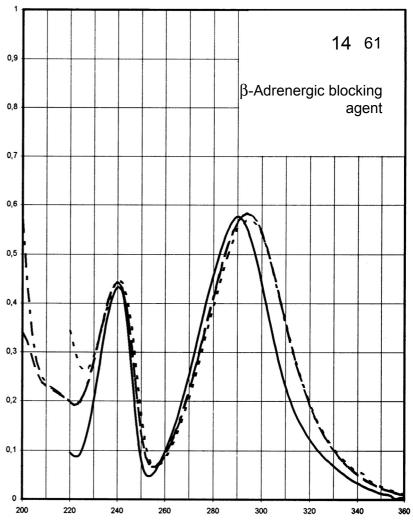

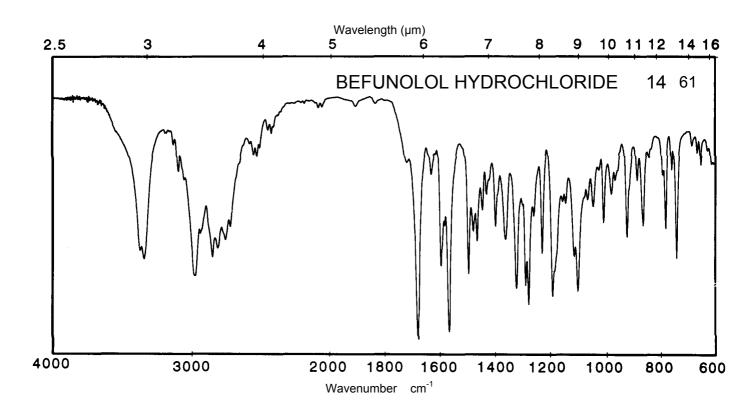





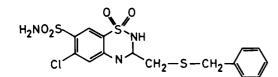

 M_r 432.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm		281 nm	296 nm
E ^{1%} _{1cm}	258		226	306
3	11160		9760	13200

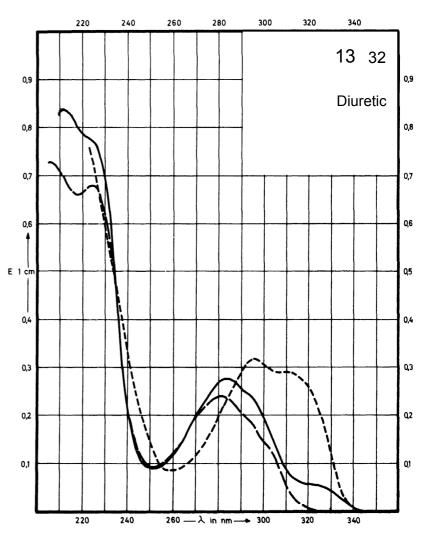

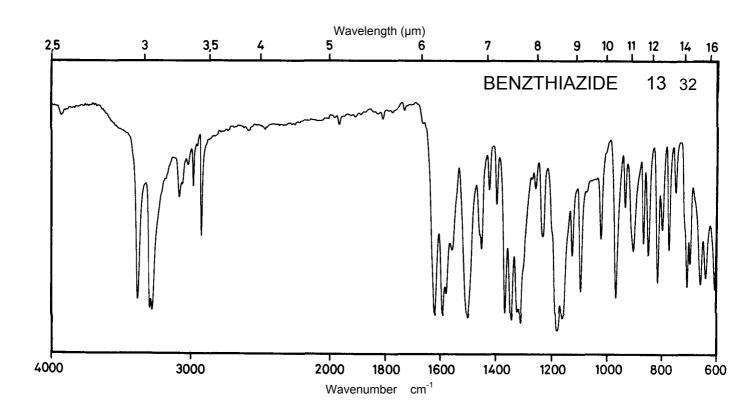



Name BEMETIZIDE

 M_r


401.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm	315 nm	315 nm	320 nm
	271 nm	272 nm	272 nm	273 nm
E ^{1%}	79	82	82	68
1cm	594	576	576	458
ε	3170	3300	3300	2700
	23900	23100	23100	18400

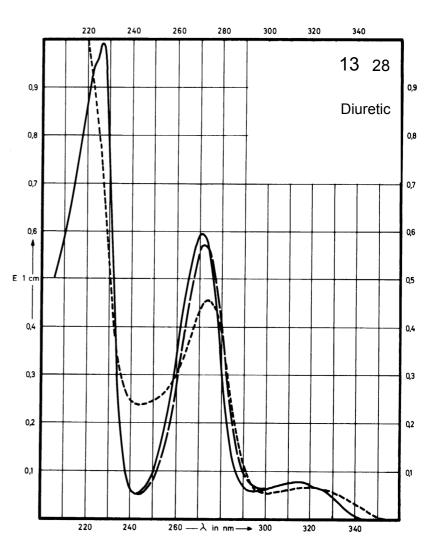

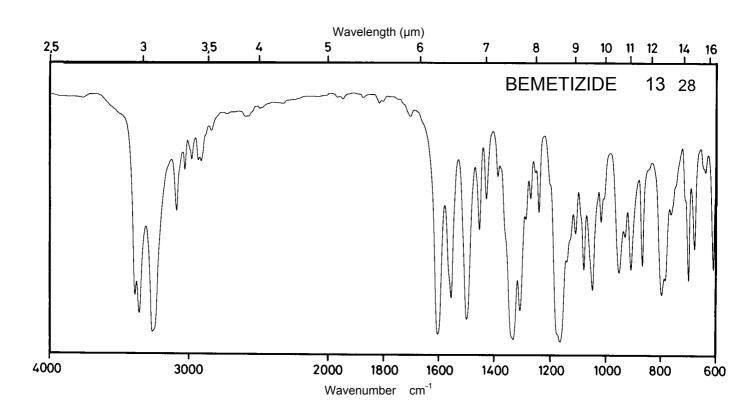




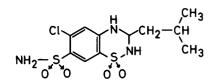


353.8

Concentration 1mg / 100 ml 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 271 nm		315 nm 272 nm	320 nm 273 nm
E ^{1%} 1cm	87 625		89 600	74 470
ε	3080 22110		3150 21230	2620 16630

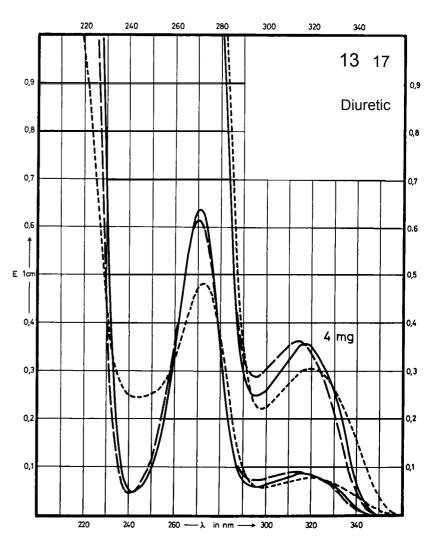

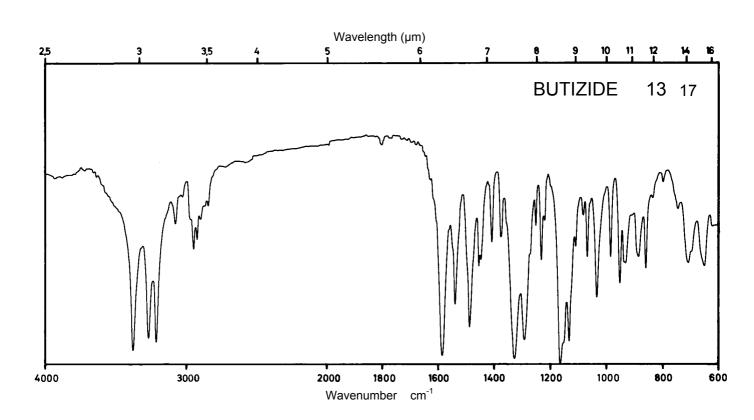


Name BENDROFLUME-THIAZIDE


$$\begin{array}{c} F_3C \\ \\ H_2N-O_2S \end{array} \begin{array}{c} H \\ S \\ O_2 \end{array} \begin{array}{c} CH_2 - CH_2 -$$

 M_{r}

421.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	324 nm 272 nm		324 nm 273 nm	330 nm 273 nm
E ^{1%} 1cm	99 560		102 535	87 432
ε	4170 23600		4300 22540	3670 18200

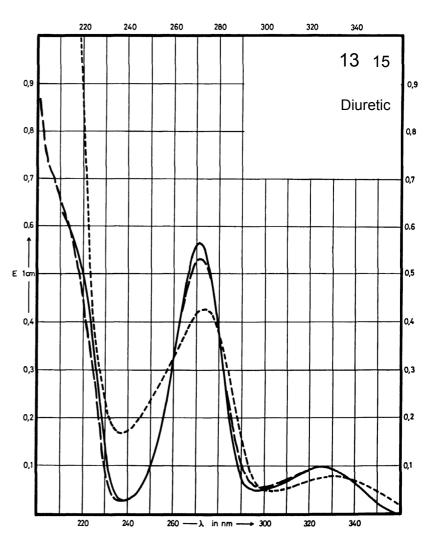

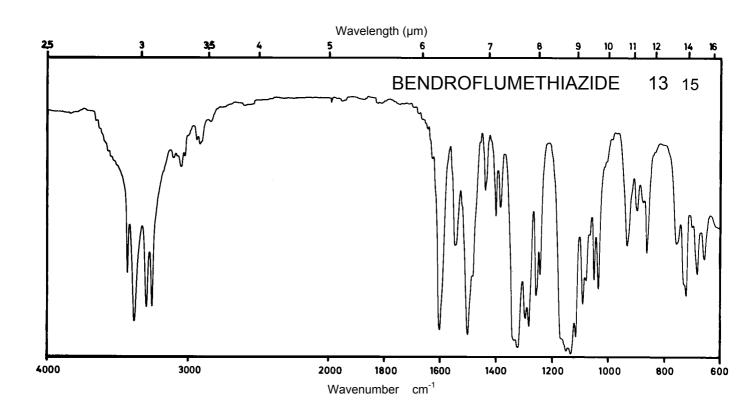

Name BUMETANIDE


H₂N-5-NH-C₄H₉


 M_{r}


364.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	335 nm 222 nm		340 nm 226 nm	317 nm
E ^{1%} 1cm	116 703		80 631	87
ε	4230 25620		2920 22990	3170



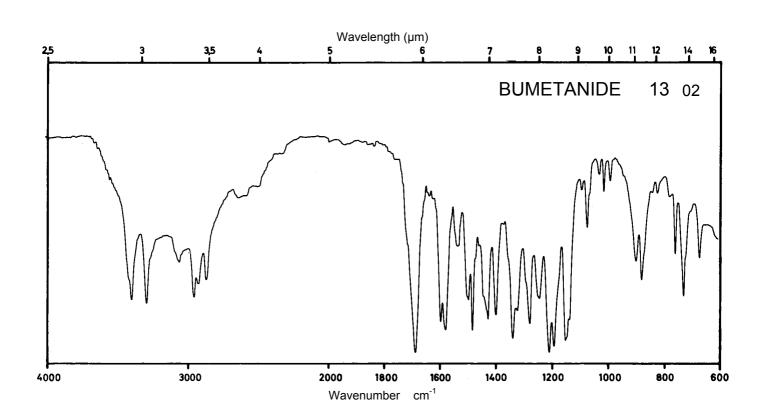

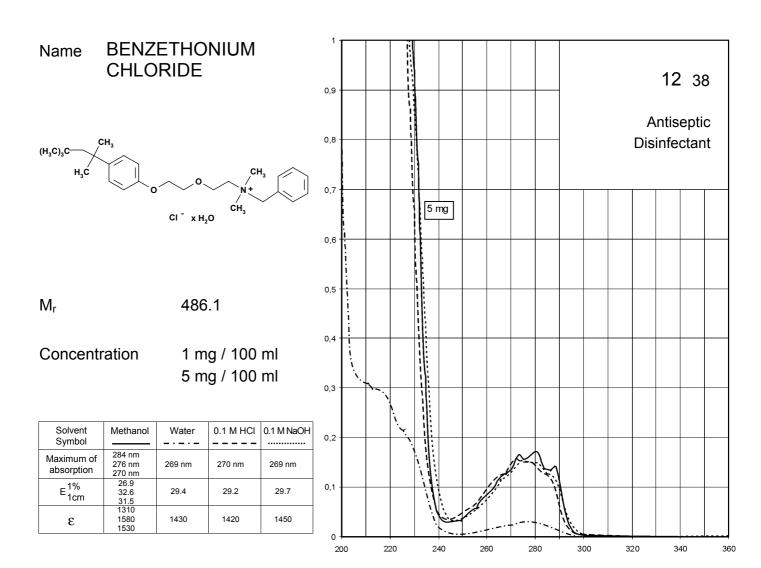

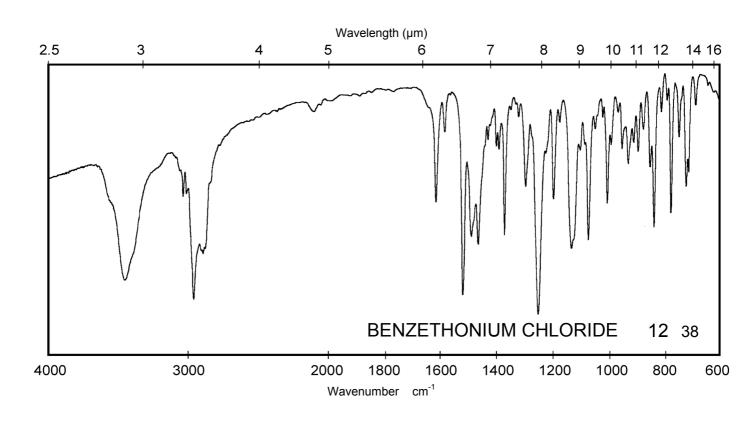



Name BUTYLPARABEN

CO-O-CH₂-CH₂-CH₂-CH₃


 M_r 194.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	256 nm		255 nm	296 nm
E ^{1%} 1cm	860		811	1245
ε	16700		15800	24200

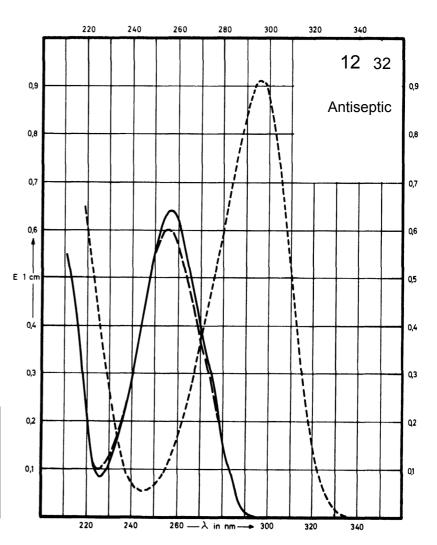

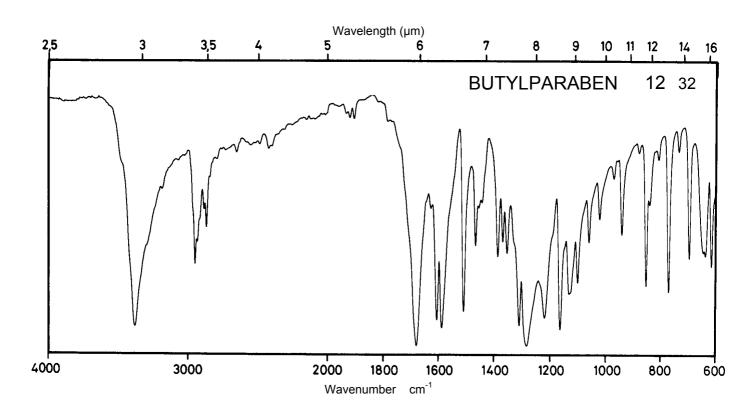





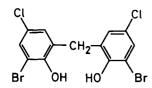

 M_{r} 426.9

Concentration 2 mg / 100 ml 4 mg / 100 ml

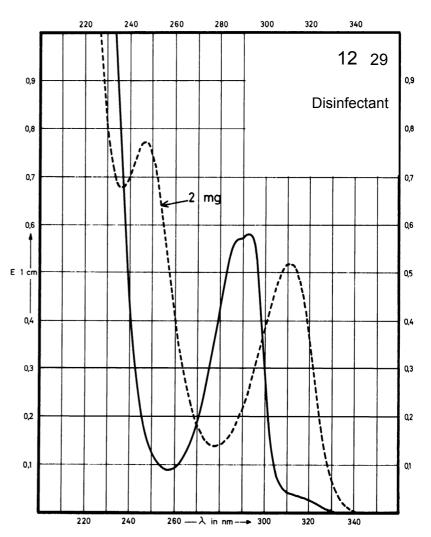


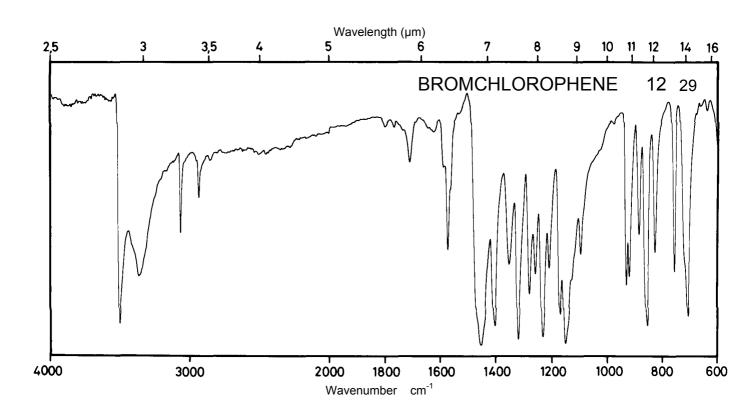


Name BENZALKONIUM CHLORIDE


 $R = C_8H_{17}$ bis $C_{18}H_{37}$

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		269 nm 263 nm 257 nm		
E ^{1%} 1cm		ca. 9.7 ca. 12.0 ca. 10.0		
3				

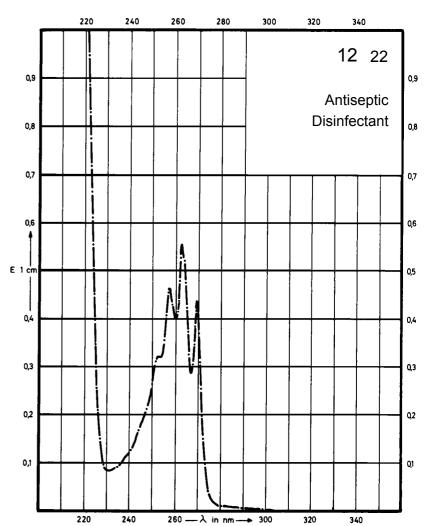

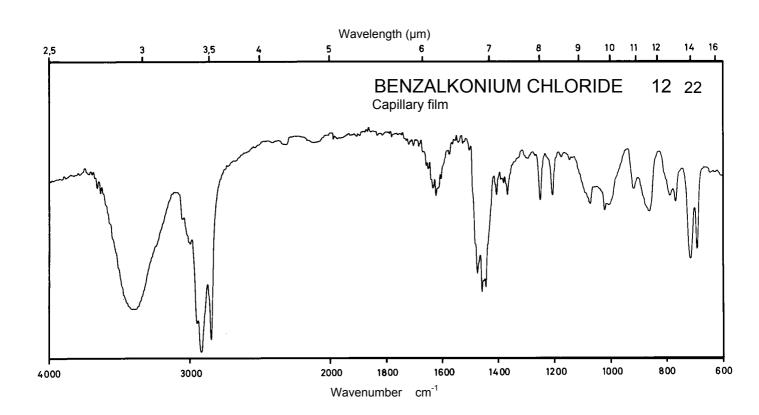





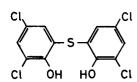

356.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	305 nm			328 nm
E ^{1%} 1cm	237			364
3	8440			12960

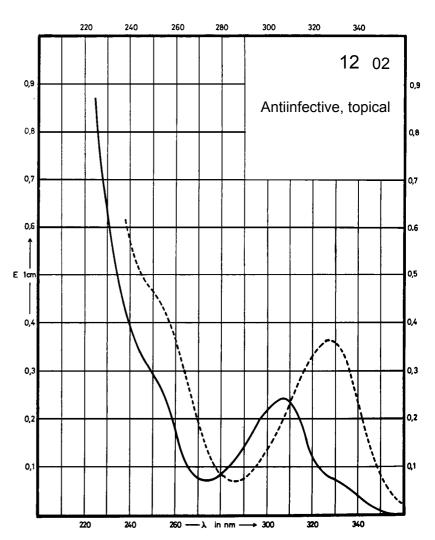

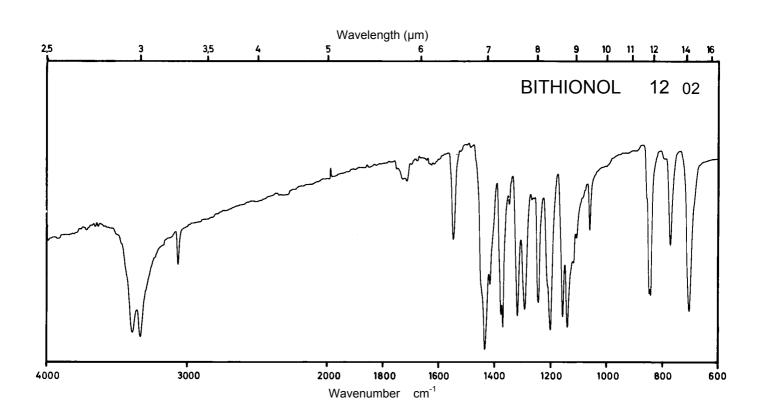



Name BENZYL ALCOHOL

CH₂OH


 M_r 108.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	252 nm 258 nm 264 nm			
E ^{1%} 1cm	14.0 17.1 13.1			
ε	151 185 142			

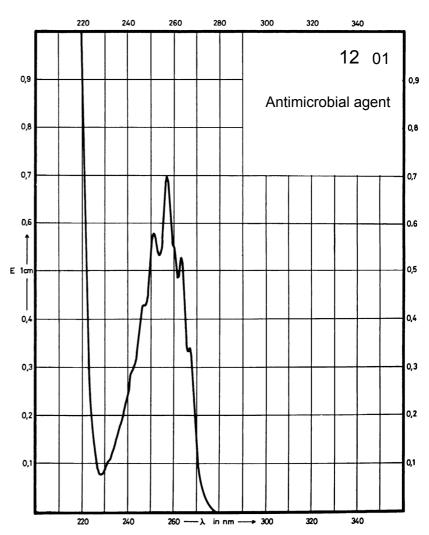

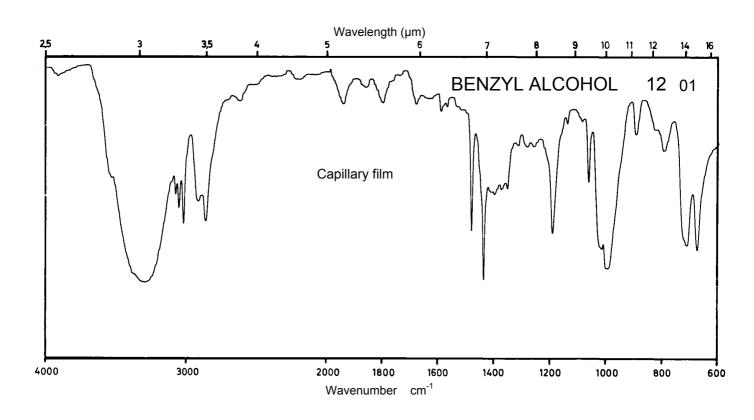



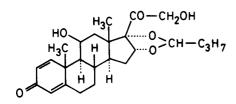


 M_r


430.5

Concentration


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm			
E ^{1%} 1cm	347			
ε	14950			

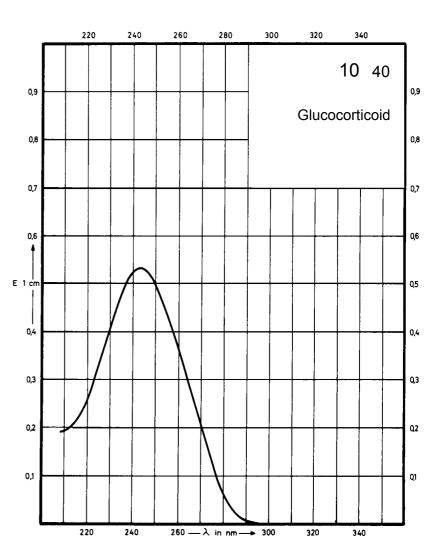

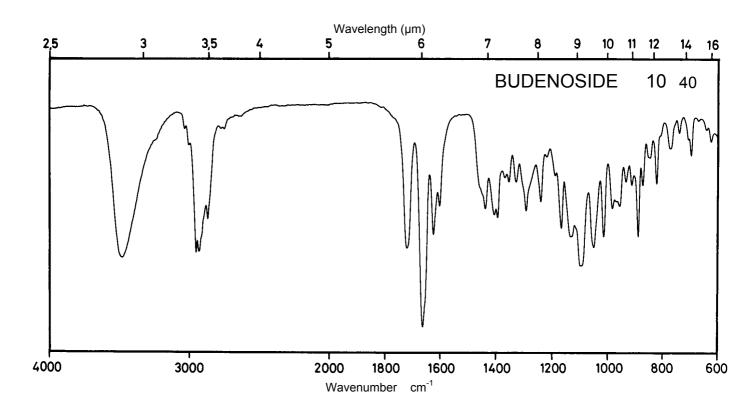




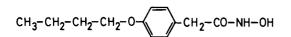


 M_r 223.3


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 284 nm 277 nm 228 nm 60 71 Maximum of 275 nm 275 nm absorption E^{1%} 1cm 59 65 1580 11500 3 1320 1440

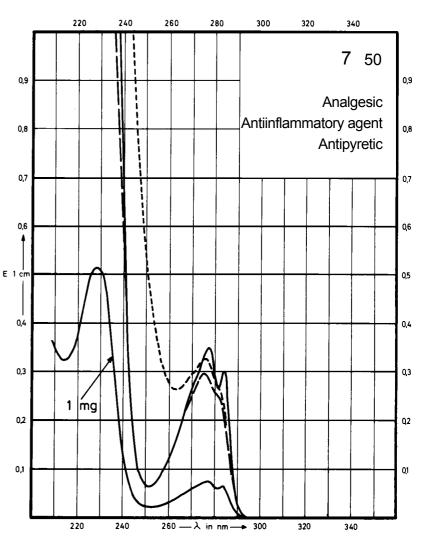

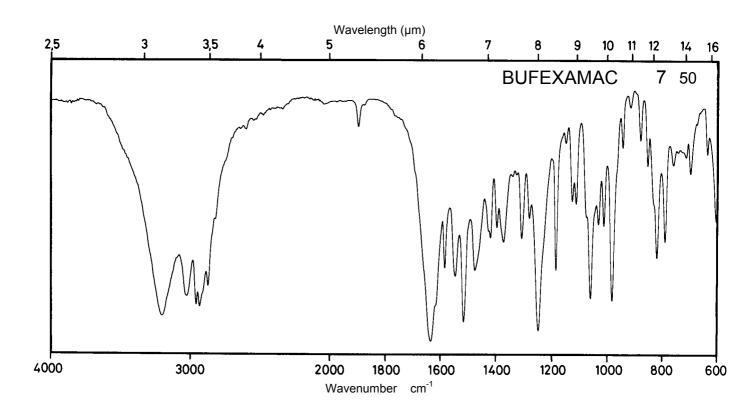


Name BETAMETHASONE 21-(DISODIUM PHOSPHATE)

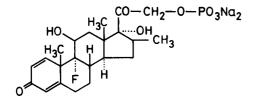


 M_r

516.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm	242 nm		
E ^{1%} 1cm	259	261		
ε	13400	13500		

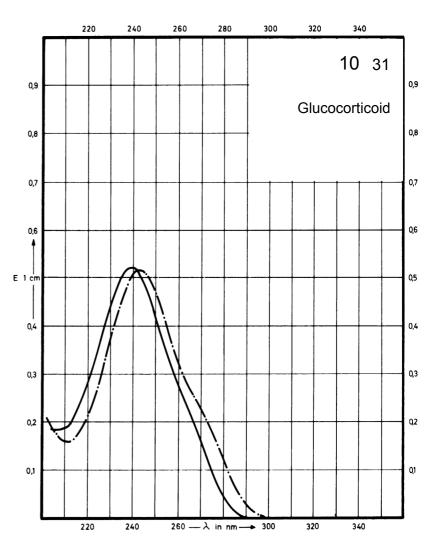

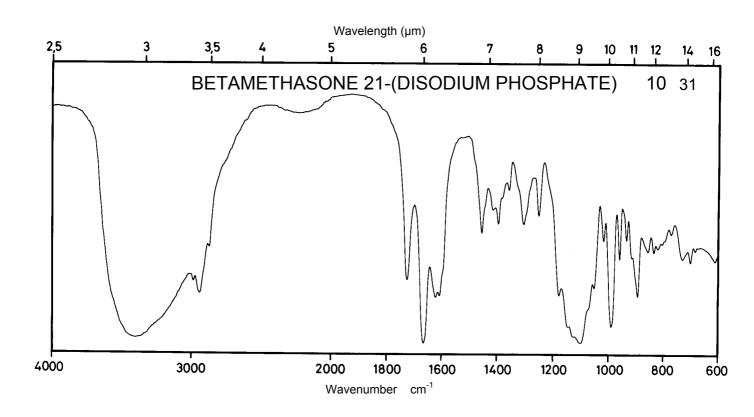


Name BECLOMETASONE 17α, 21-DIPRO-PIONATE


$$CH_2-O-CO-CH_2-CH_3$$
 CO
 CO
 CH_3
 CH_3

 M_r


521.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm			
E ^{1%} 1cm	296			
ε	15400			

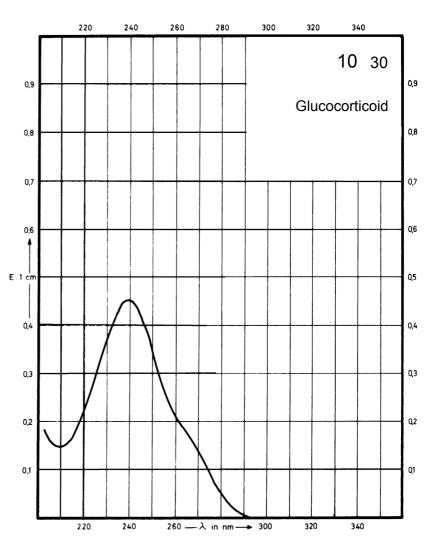

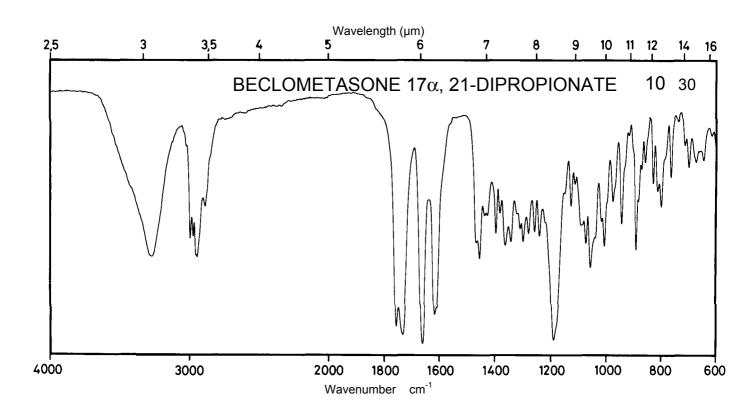





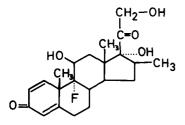


392.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm			
E ^{1%} 1cm	387			
ε	15190			

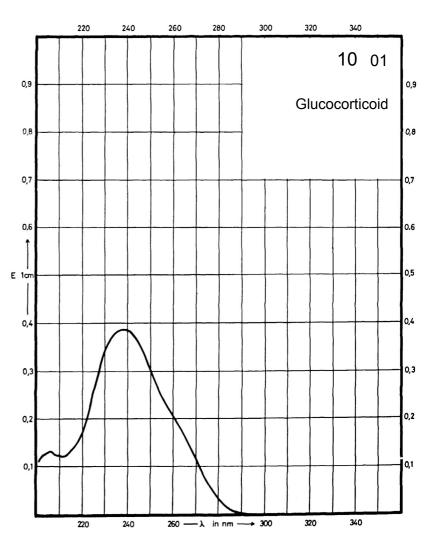

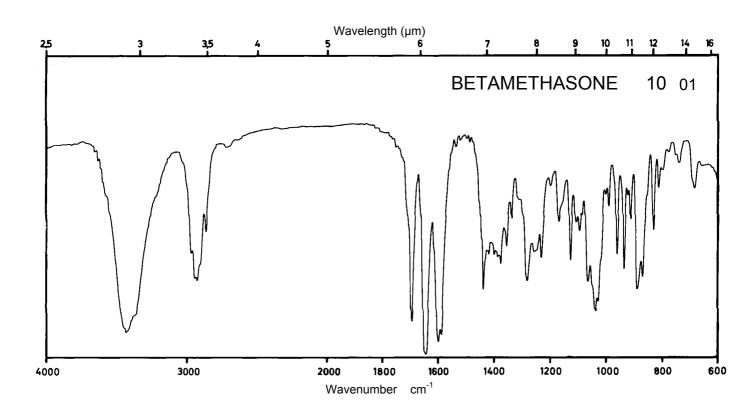





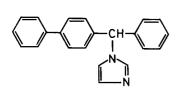


310.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	254 nm		255 nm	
E ^{1%} 1cm	778		772	
ε	24200		24000	

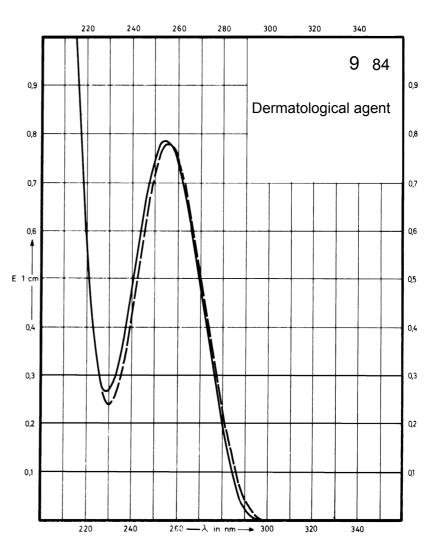

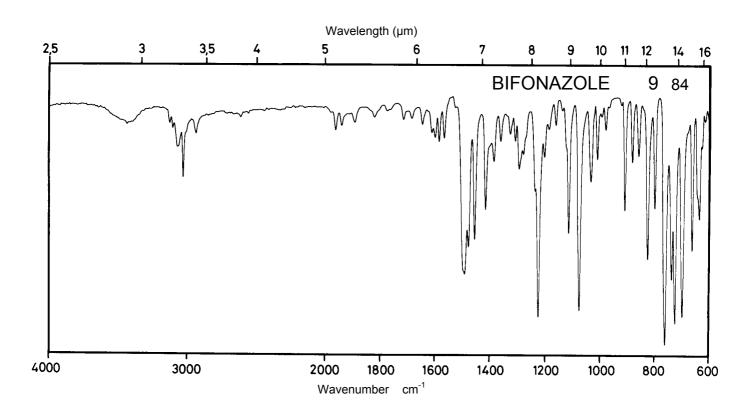





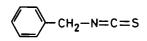



149.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	246 nm			
E ^{1%} 1cm	104			
ε	1550			

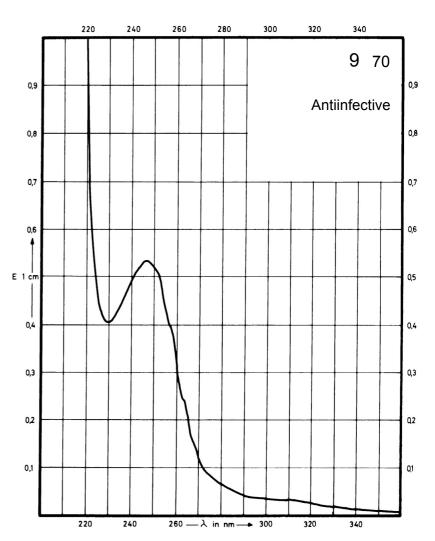

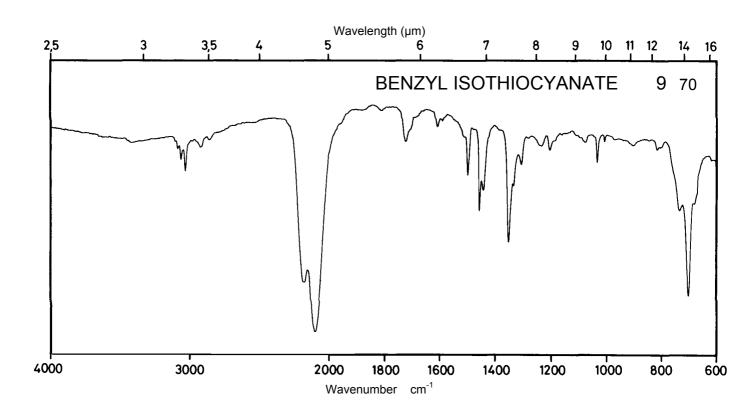

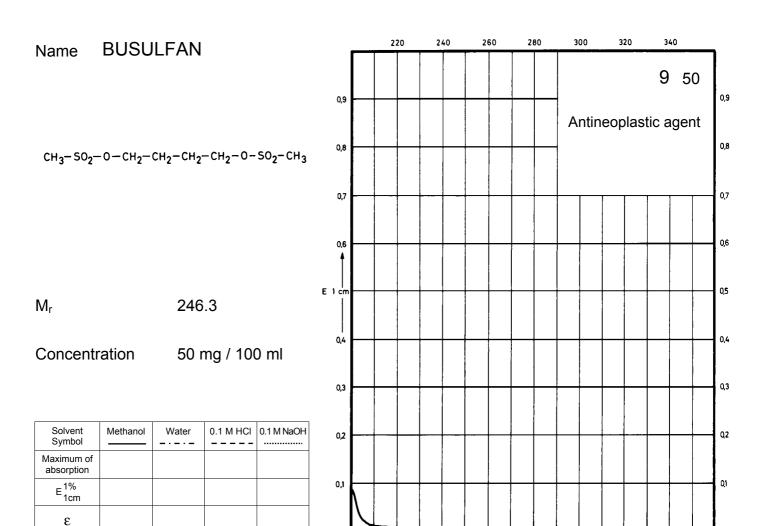

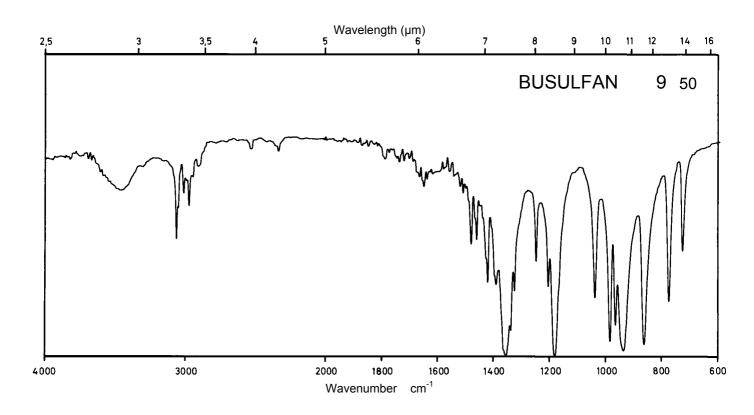


220


Name BENSULDAZIC ACID SODIUM SALT

 M_{r}

322.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm 248 nm		285 nm 251 nm	290 nm 258 nm
E ^{1%} 1cm	354 267		359 288	405 440
3	11410 8610		11570 9290	13060 14190

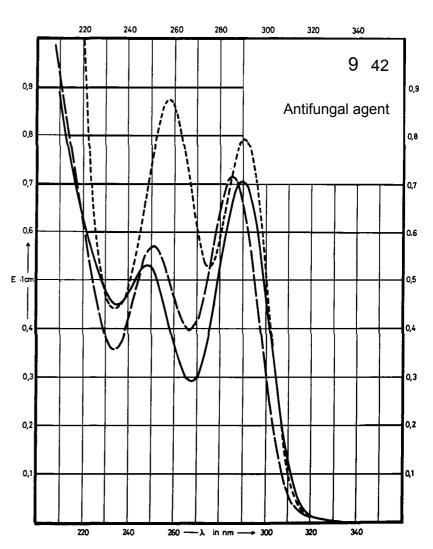

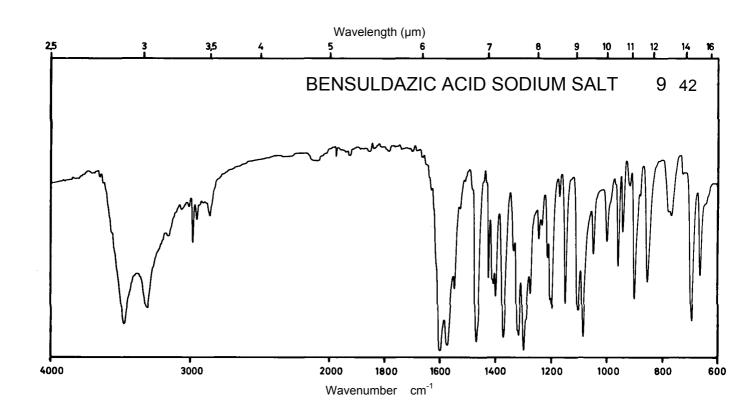


Name BUCLOSAMIDE


$$CI - CO-NH-CH_2-CH_2-CH_2-CH_3$$


 M_{r}


227.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	301 nm 244 nm		297 nm 245 nm	324 nm 246 nm
E ^{1%} 1cm	233 521		204 532	322 391
3	5300 11860		4640 12110	7330 8900

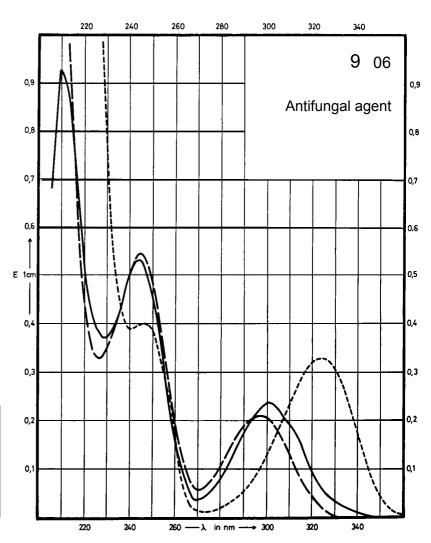


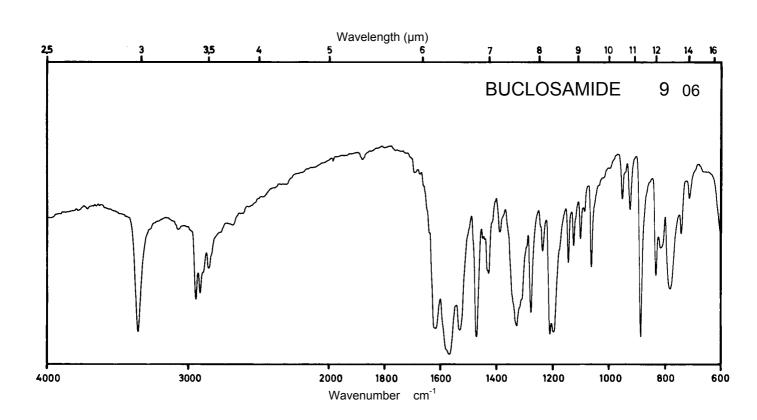

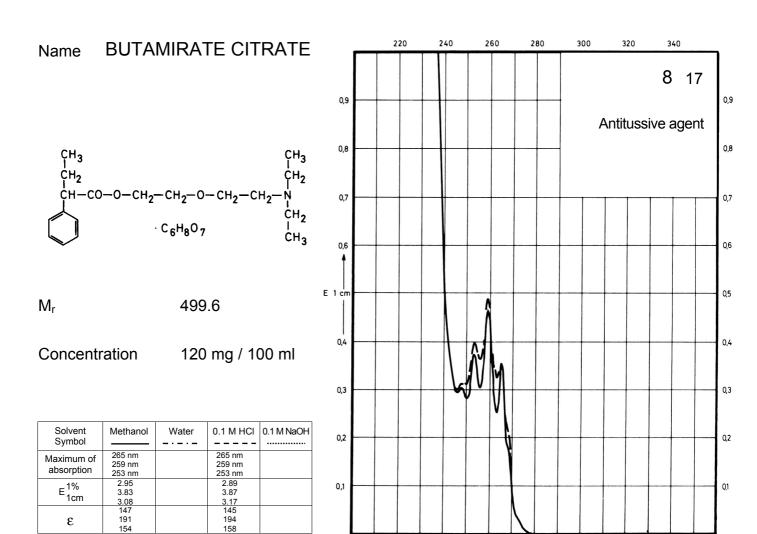

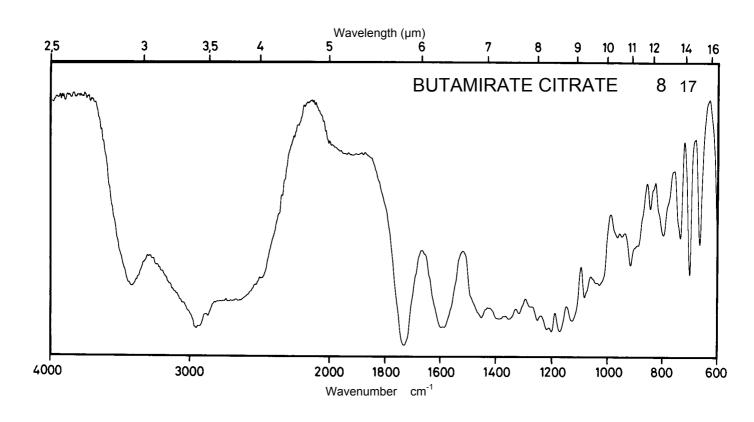

Name BROMHEXINE HYDROCHLORIDE

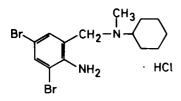



 M_{r}

412.6


Concentration 2 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	318 nm 250 nm		312 nm 247 nm	
E ^{1%} 1cm	83 266		75 233	
ε	3420 10990		3090 9600	

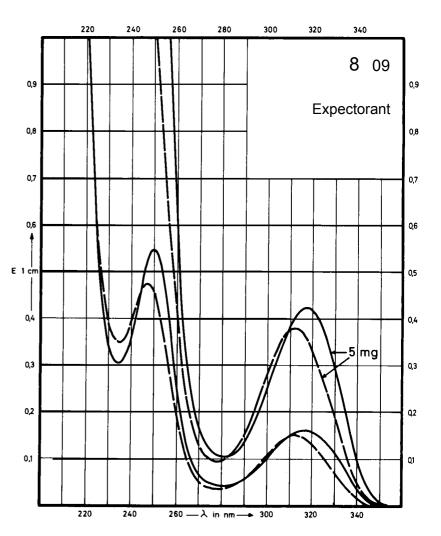

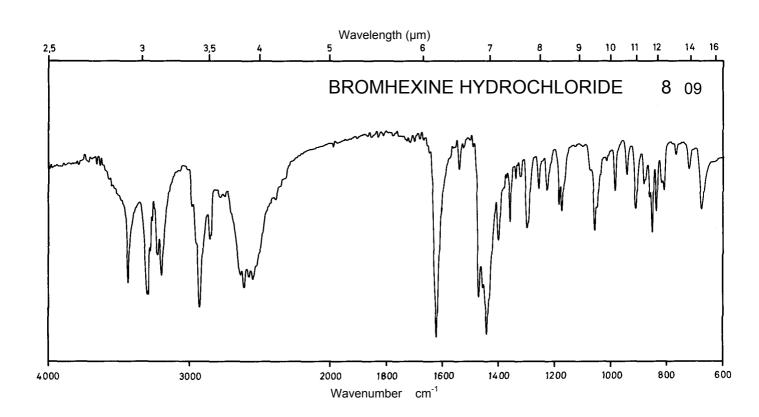




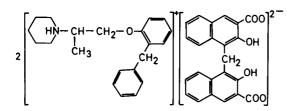


 M_r 1007.3

Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	368 nm 278 nm 237 nm			
E ^{1%} 1cm	66 115 1020			
ε	6600 11630 103000			

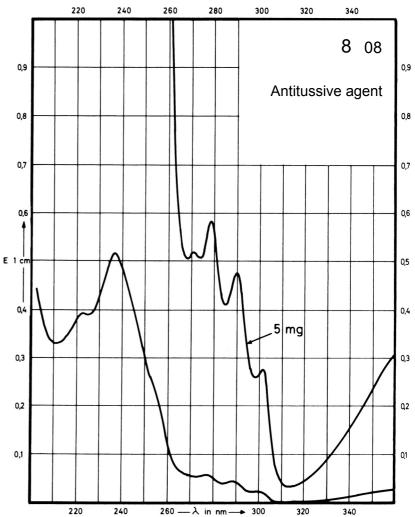

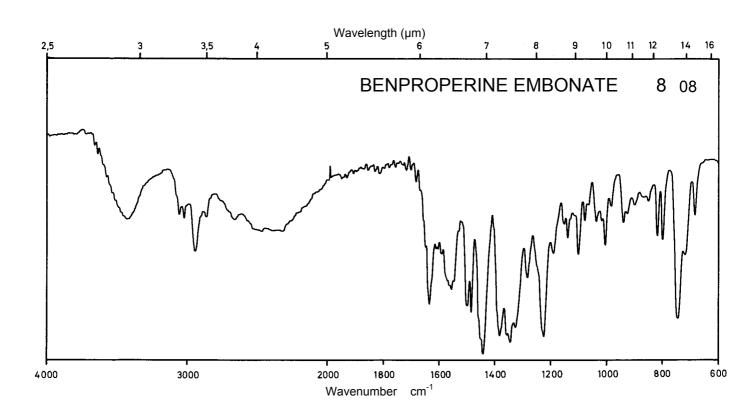



Name BENORILATE

 M_{r}


313.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm		239 nm	Decom- position observed
E ^{1%} _{1cm}	747		743	
ε	23400		23300	

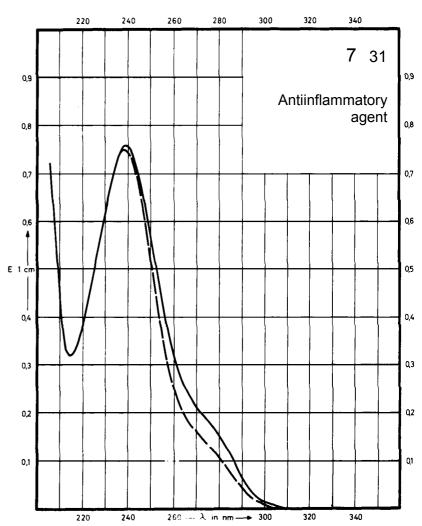

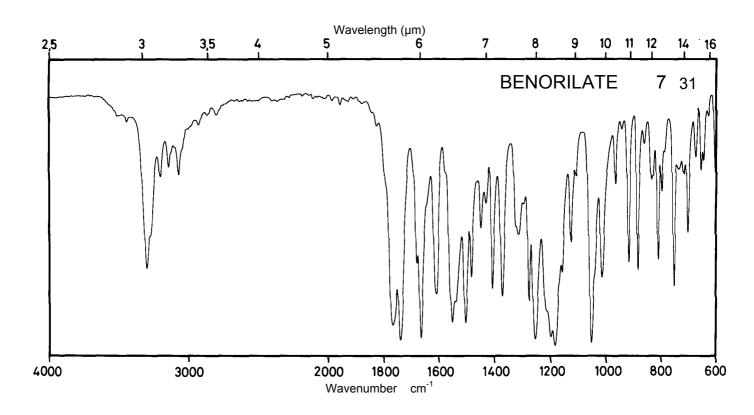




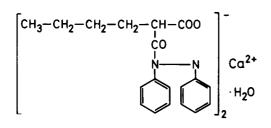


 M_{r} 690.9


	Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
	Maximum of absorption	237 nm	234 nm	233 nm	234 nm
	E ^{1%} *	509	447	454	449
	ε	35100	30900	31400	31000

^{*} Calculated on dried substance

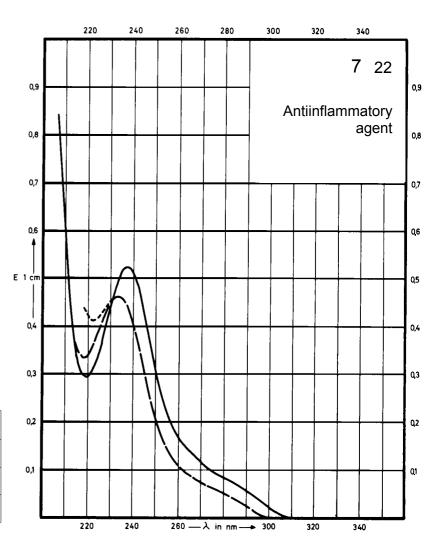

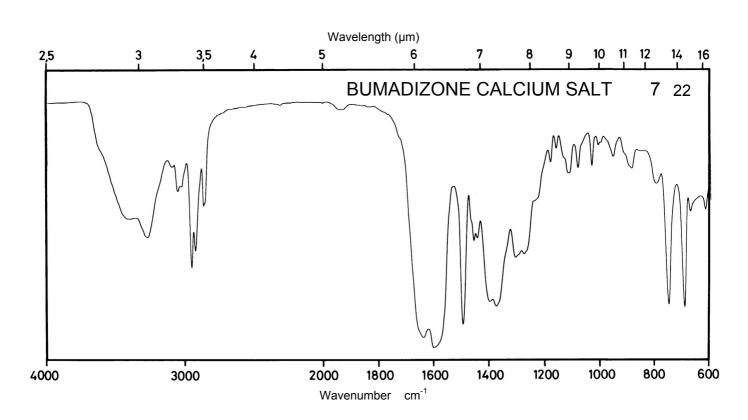





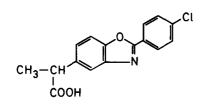


301.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	308 nm 243 nm			309 nm 244 nm
E ^{1%} 1cm	945 302			900 295
3	28500 9100			27100 8900

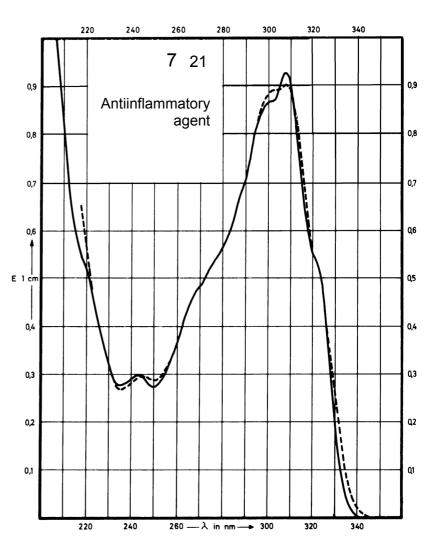

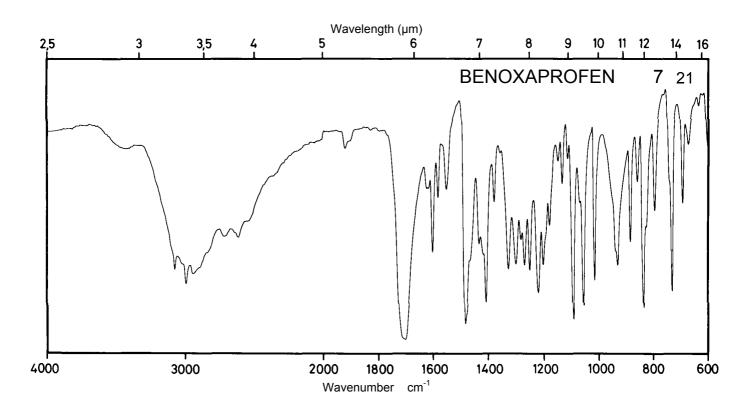



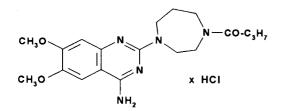



409.9

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	343 nm 330 nm 246 nm	331 nm 245 nm	330 nm 245 nm	347 nm 275 nm 250 nm
E ^{1%} 1cm	242 259 1235	249 1190	251 1182	137 370 1295
3	9050 9700 46200	9300 44400	9400 44200	5100 13800 48400

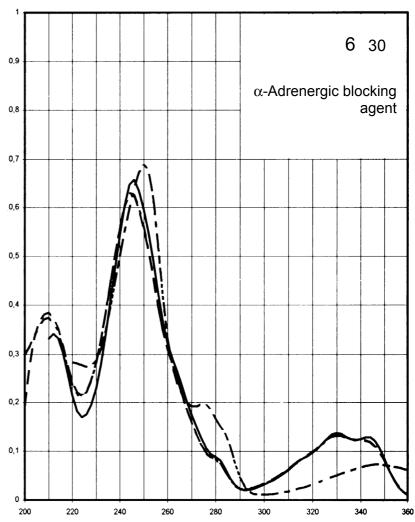

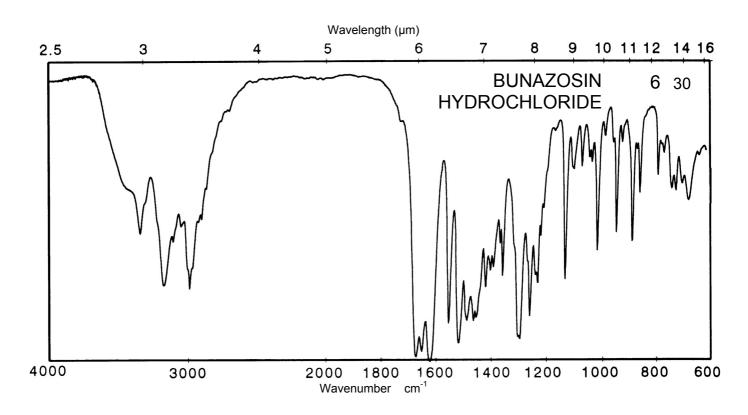





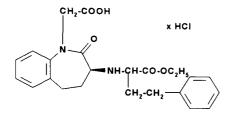

461.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm	240 nm	237 nm	241 nm
E ^{1%} 1cm	219	218	213	222
3	10100	10000	9800	10200

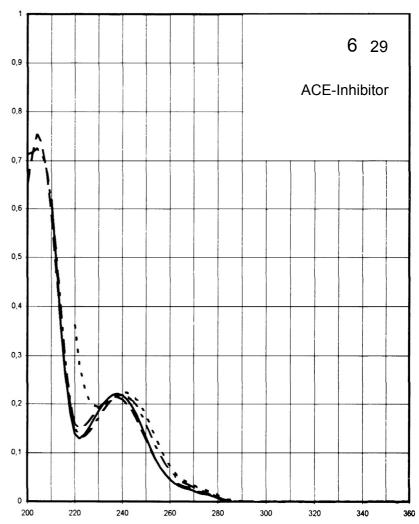

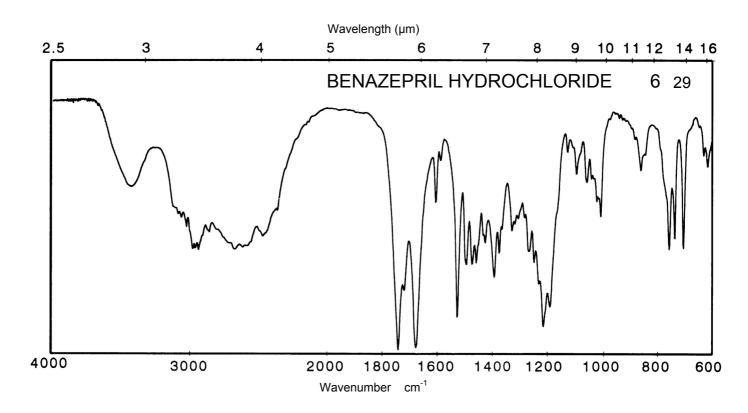



Name BUFORMIN HYDROCHLORIDE

 M_{r}

193.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		232 nm		230 nm
E ^{1%} _{1cm}		774		746
3		14990		14450

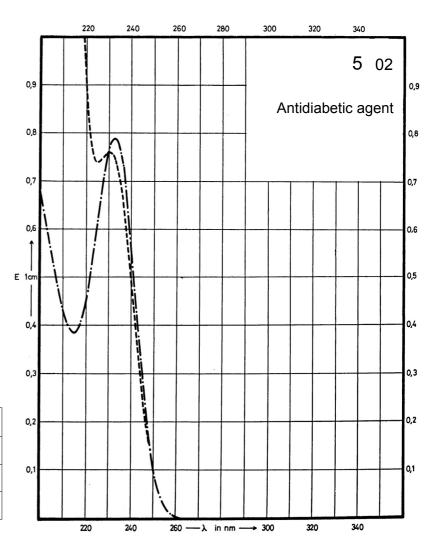

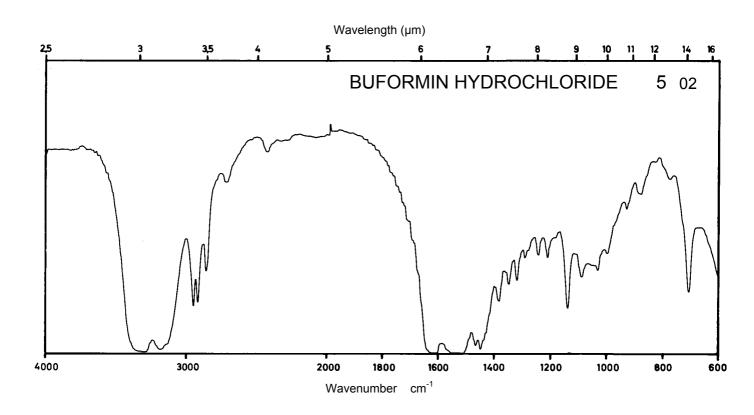




 $C_{16}H_{20}N_2$. $(C_{16}H_{18}N_2O_4S)_2$


 M_r

909.1


Concentration 1 mg / 100 ml 53 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm			
E ^{1%} 1cm	8.9			
3	810			

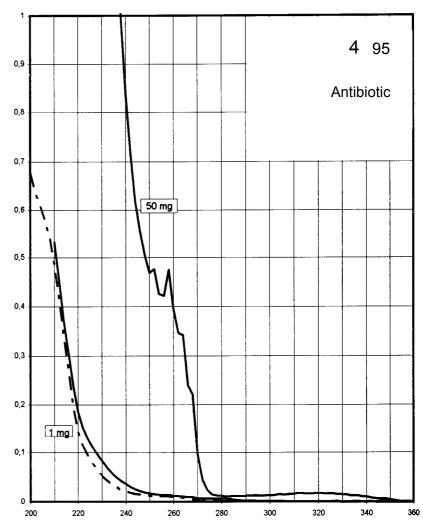

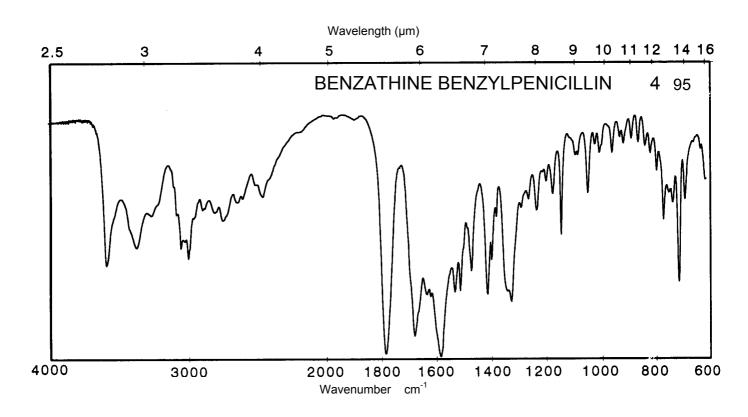




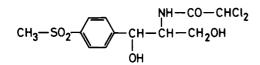


 M_r 356.2


Concentration 4 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 265 nm		273 nm 266 nm	
E ^{1%} 1cm	20.0 22.5		22.7 25.7	
ε	710 800		810 920	

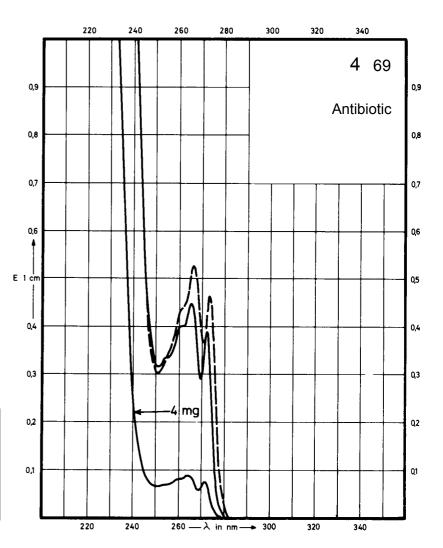

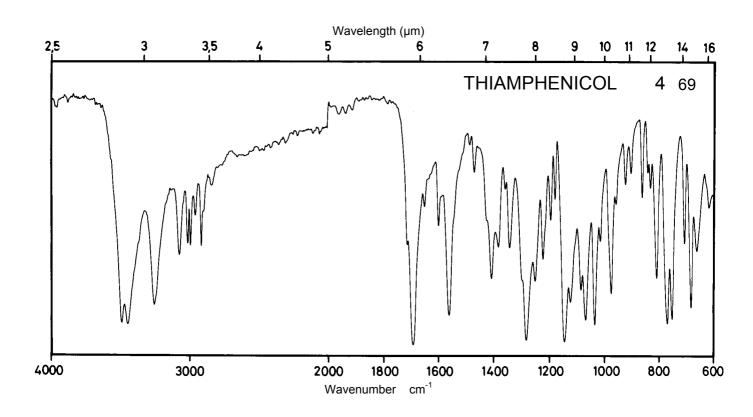


Name BLEOMYCIN SULFATE

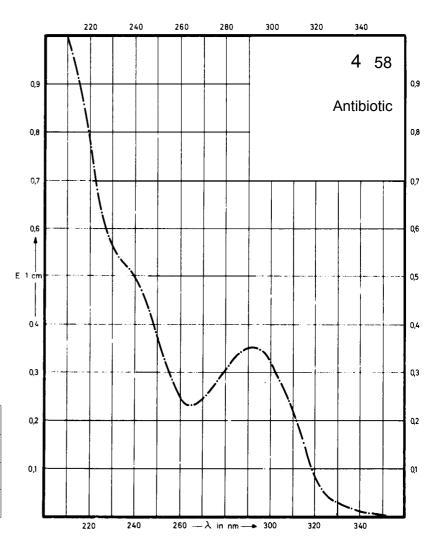


Concentration

3.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			291 nm	
E ^{1%} 1cm			ca. 100	
ε				

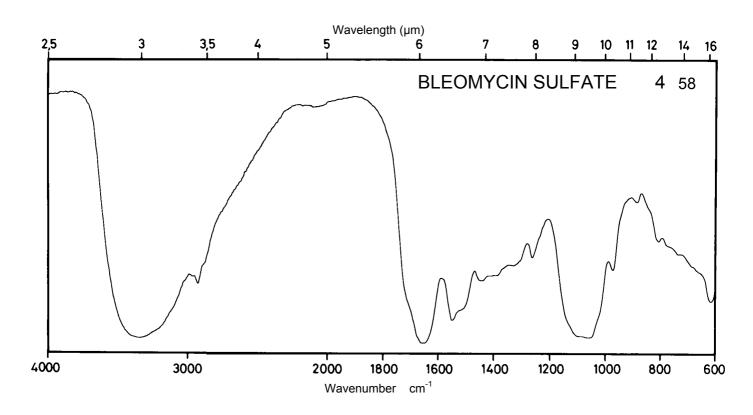





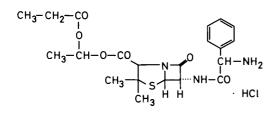


 M_r 502.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm 263 nm 258 nm	269 nm 263 nm 258 nm		
E ^{1%} 1cm	7.6 10.3 11.7	8.9 12.1 13.1		
ε	380 520 590	445 605 655		

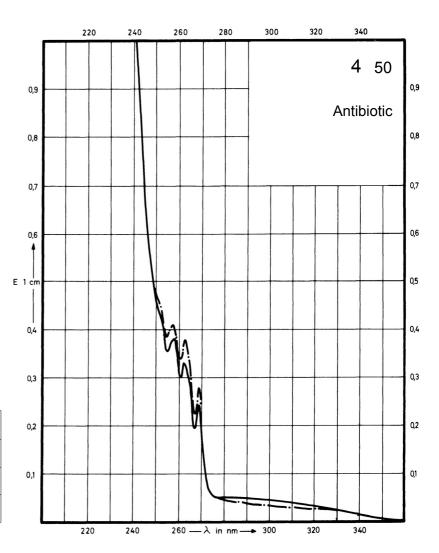

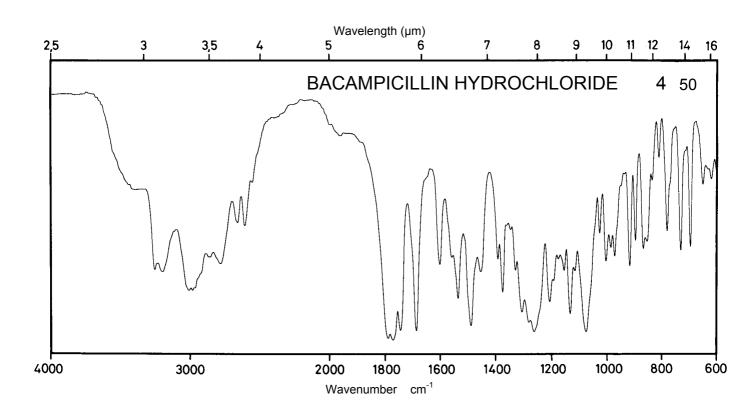


Name BACITRACIN ZINC SALT

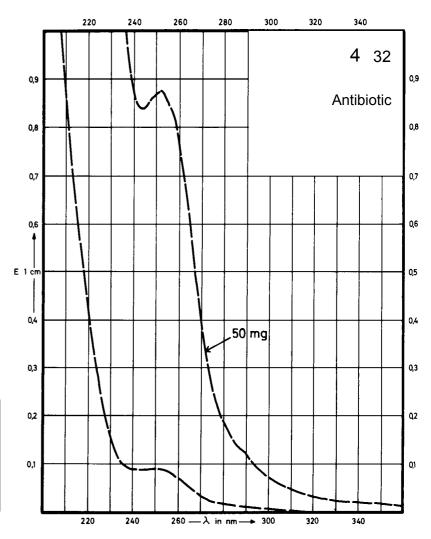


Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of absorption			252 nm	
E ^{1%} 1cm			18	
ε				

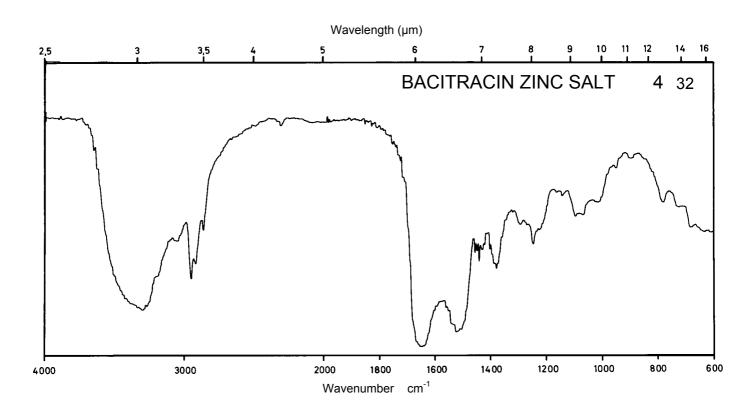





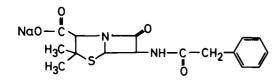


 M_r 356.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 257 nm	264 nm 257 nm		
E ^{1%} 1cm	4.6 6.8	4.7 7.1		
ε	165 242	168 253		

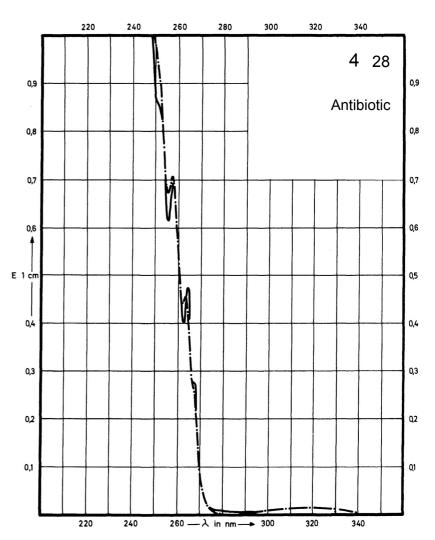

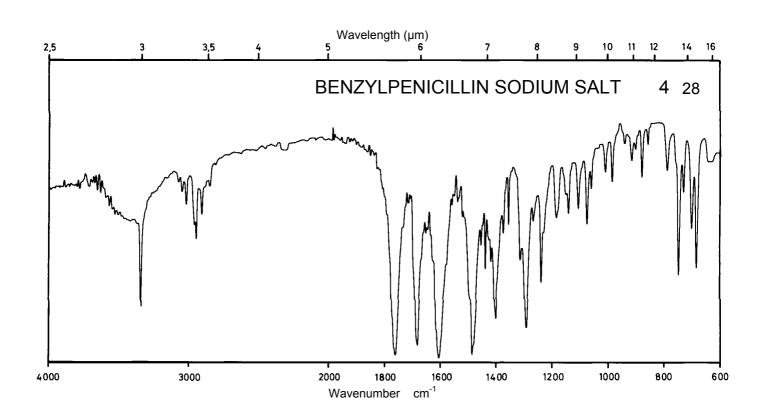




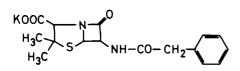


 M_r


372.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 258 nm	265 nm 258 nm		
E ^{1%} 1cm	4.6 6.6	4.4 6.8		
ε	170 245	160 250		

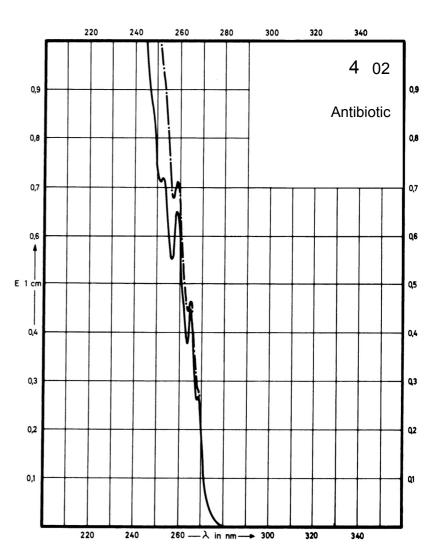

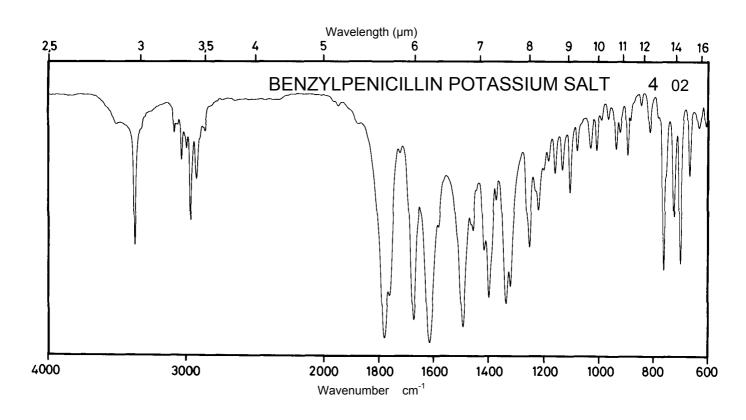




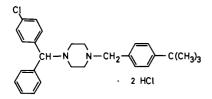


 M_{r} 506.0


Concentration 1 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm 220 nm		224 nm	
E ^{1%} 1cm	17.2 431		434	
ε	870 21800		21950	

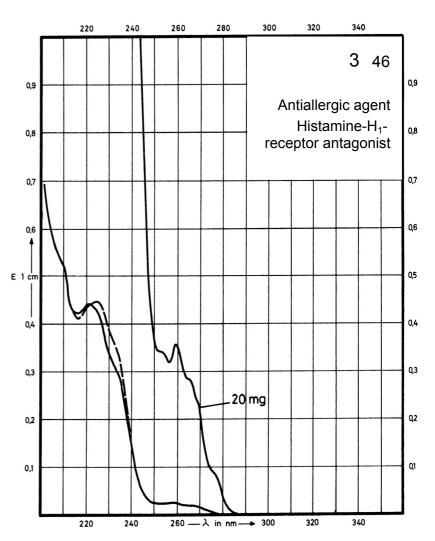

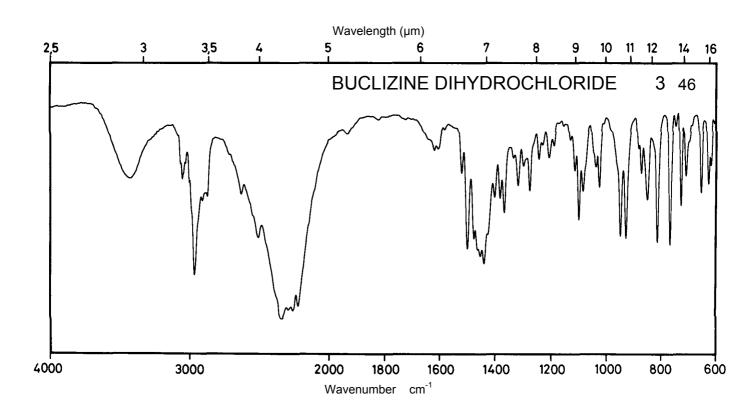



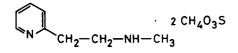



 M_r

Concentration 3.0 mg / 100 ml 4.6 mg / 100 ml


328.4

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 260 nm 261 nm 261 nm absorption E^{1%} 1cm 104 223 117 3 3400 7320 3850

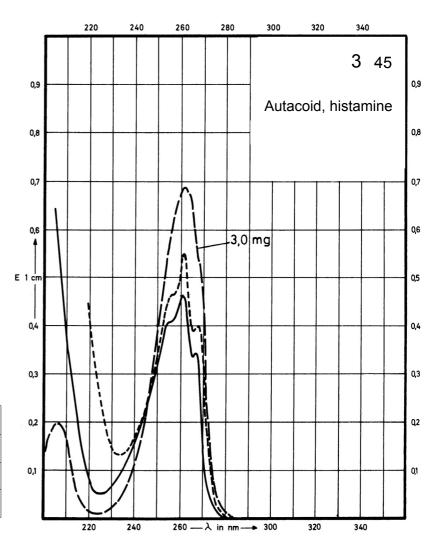

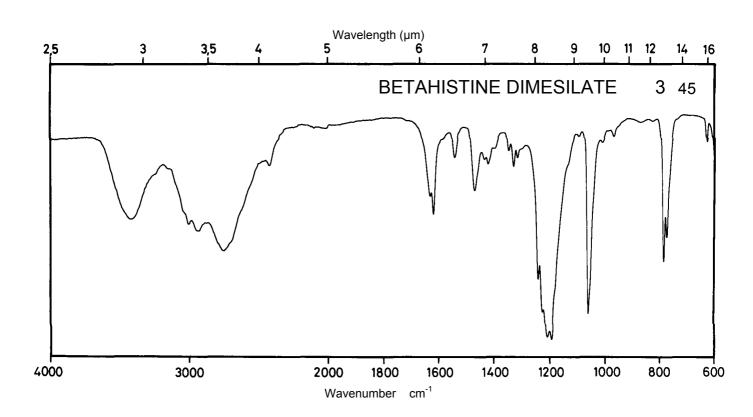


Name BETAHISTINE DIHYDROCHLORIDE

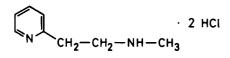


 M_r 209.1


Concentration 1.9 mg / 100 ml 2.9 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 260 nm 261 nm 261 nm absorption E^{1%} 1cm 163 350 184 ε 3400 7320 3850

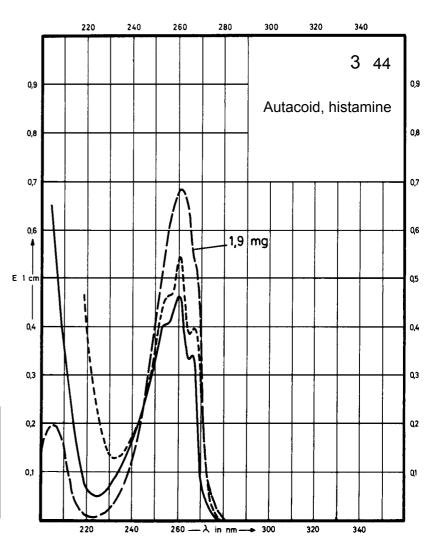

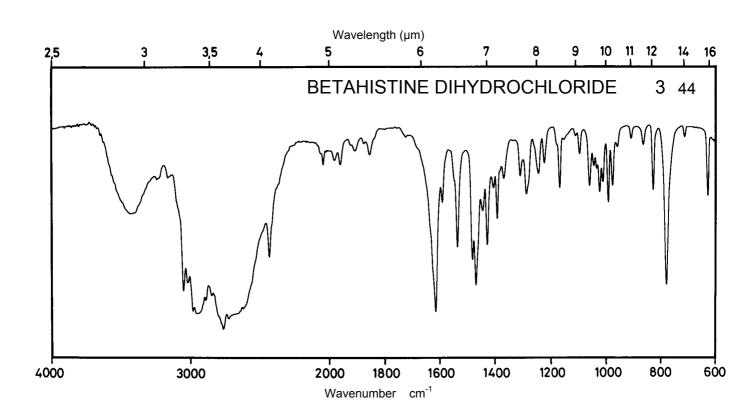





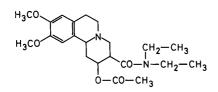


404.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm		282 nm	283 nm
E ^{1%} 1cm	94		89	93
ε	3800		3600	3770

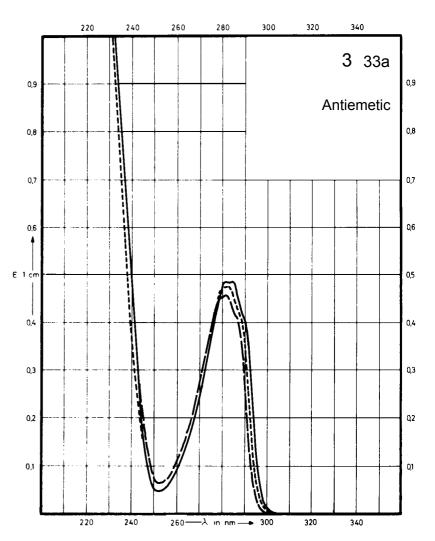

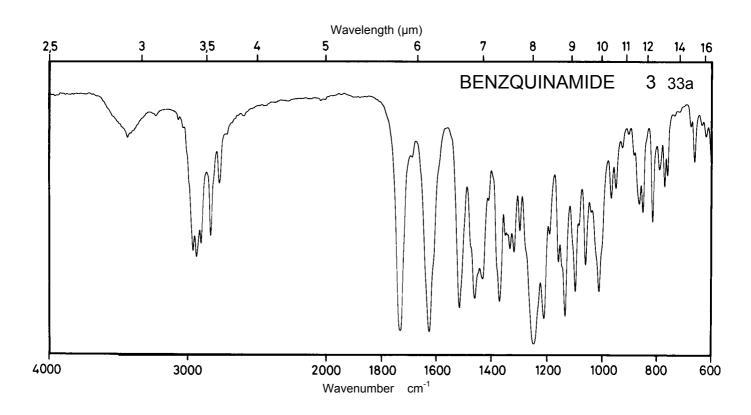





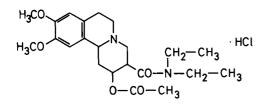

 M_{r} 441.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm		282 nm	283 nm
E ^{1%} 1cm	84		81	85
3	3700		3600	3700

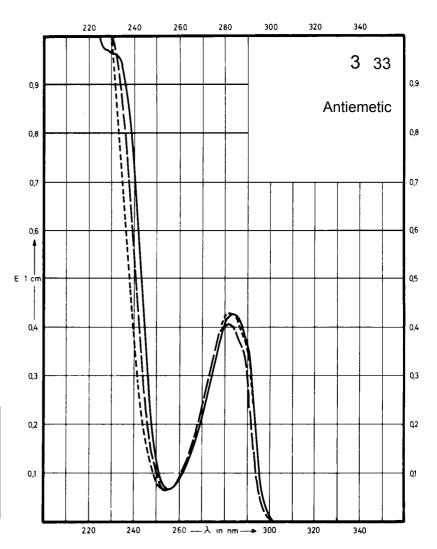

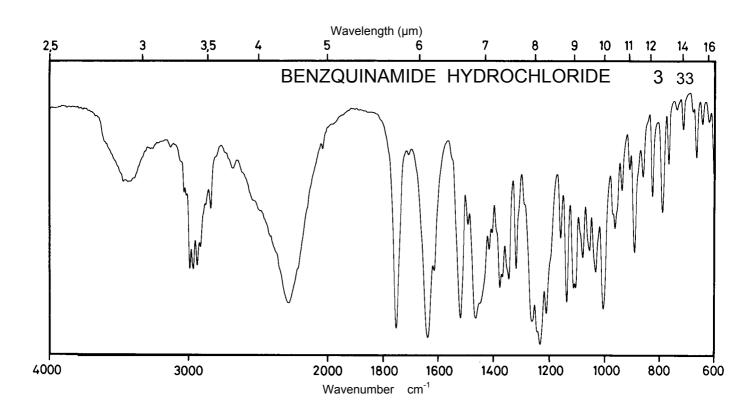





$$\begin{array}{c} \operatorname{Br} & \operatorname{CH-CH_2-CH_2-N} & \operatorname{CH_3} \\ \operatorname{CH_3} & \operatorname{C_4H_4O_4} \end{array}$$


435.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm	262 nm	264 nm	262 nm
E ^{1%} 1cm	131	137	198	132
3	5680	5980	8640	5730

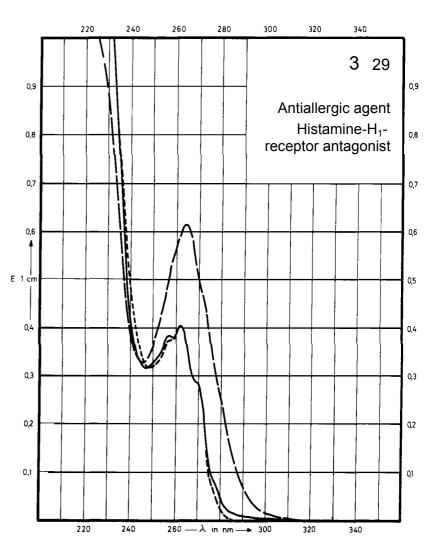

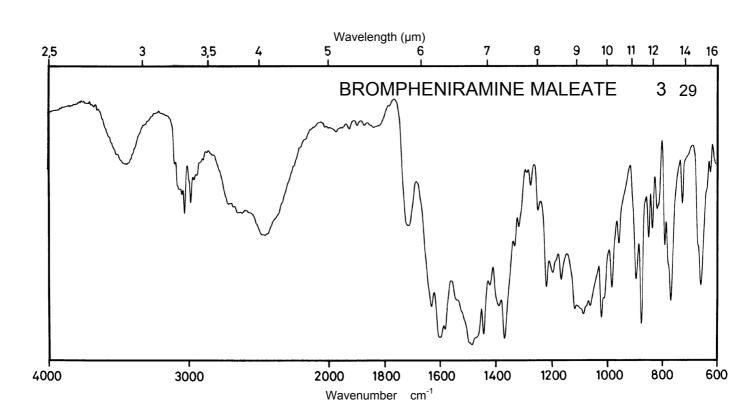




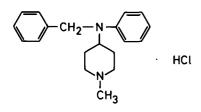


316.9


Concentration 2 mg / 100 ml 10 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	298 nm 251 nm		257 nm	
E ^{1%} 1cm	62.5 439		15.5	
3	1980 13920		490	

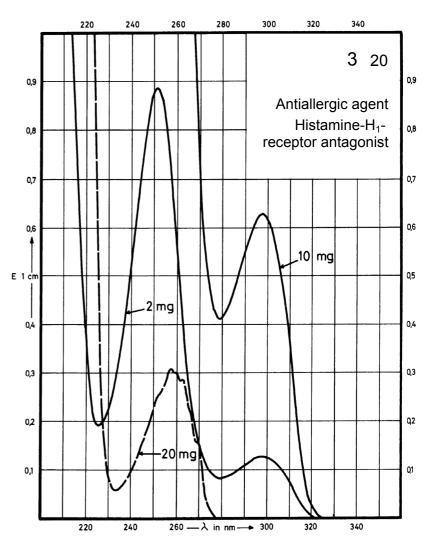

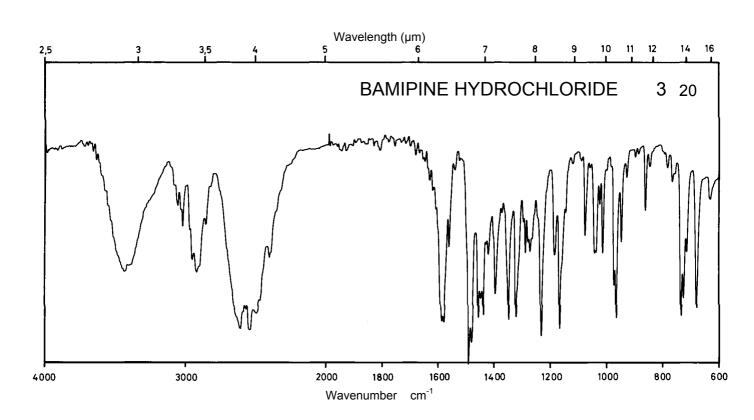

Name BROMAZINE HYDROCHLORIDE

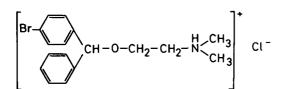



 M_r

370.7

Concentration 1 mg / 100 ml 40 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm 227 nm		258 nm 229 nm	
E ^{1%} 1cm	14.7 420		16.9 425	
ε	545 15570		625 15750	

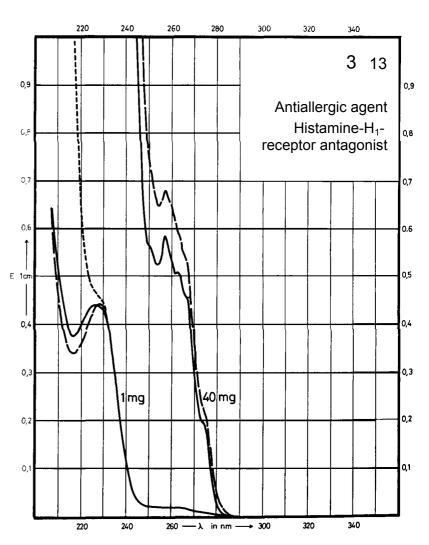

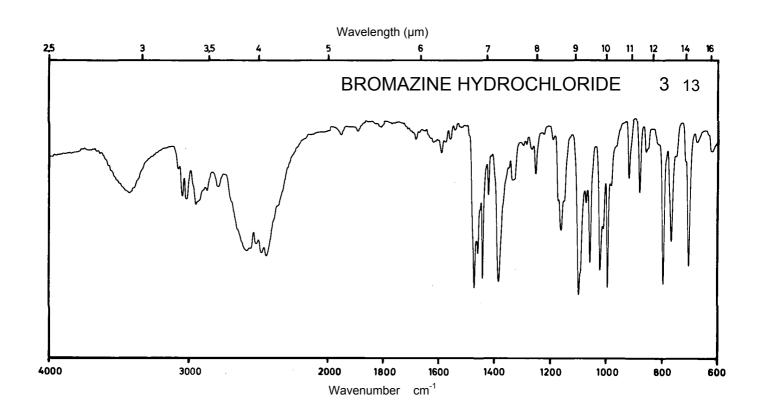

Name BROTIANIDE


$$CI \longrightarrow CS-NH \longrightarrow Br$$


$$O-CO-CH_3$$

 M_{r}


463.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	318 nm			390 nm 278 nm
E 1% 1cm	266			288 230
ε	12330			7560 10660

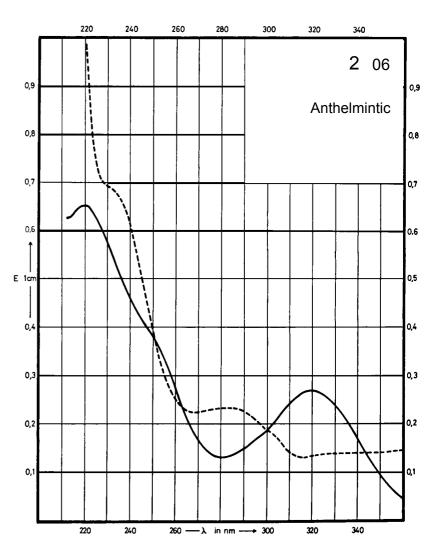

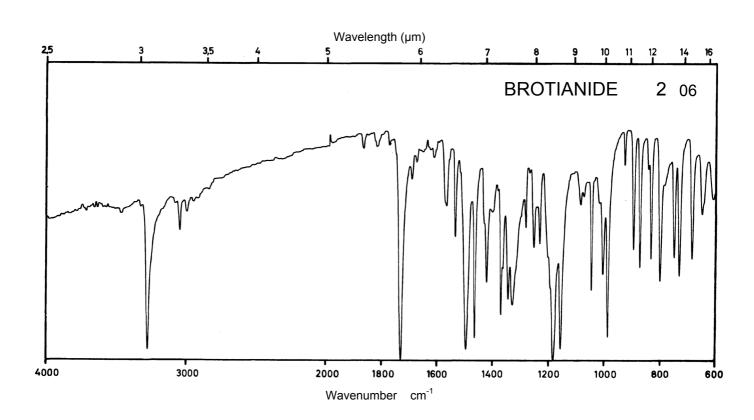




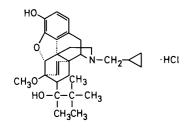


 M_r 504.1


Concentration 4 mg / 100 ml 12 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm	286 nm	286 nm	299 nm
E ^{1%} 1cm	35	31	31	51
ε	1760	1560	1560	2580

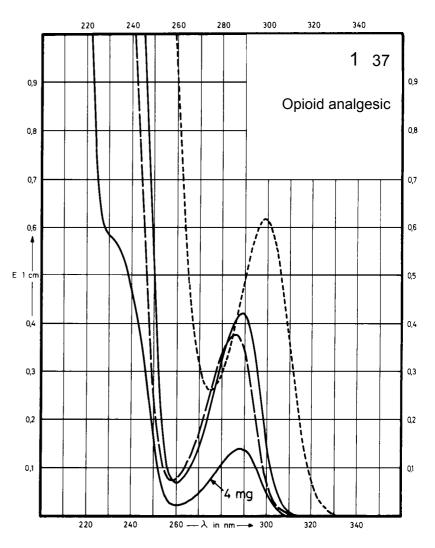

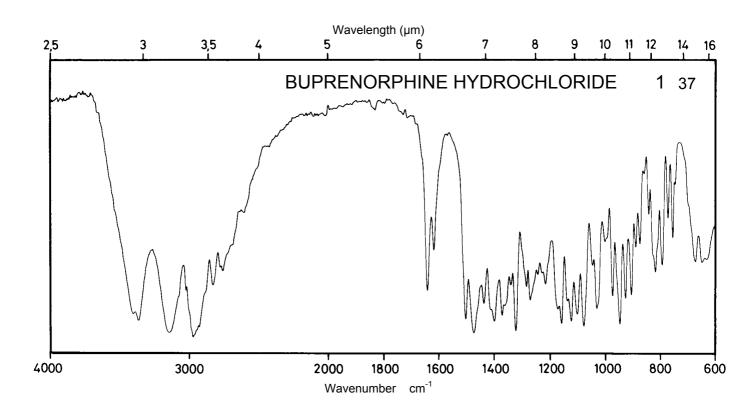



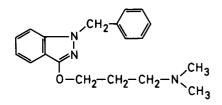



309.4

Concentration


1.9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	307 nm		306 nm	306 nm
E ^{1%} 1cm	167		167	161
ε	5170		5170	4980

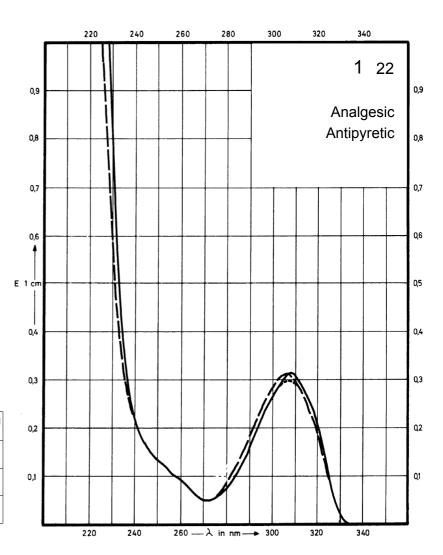

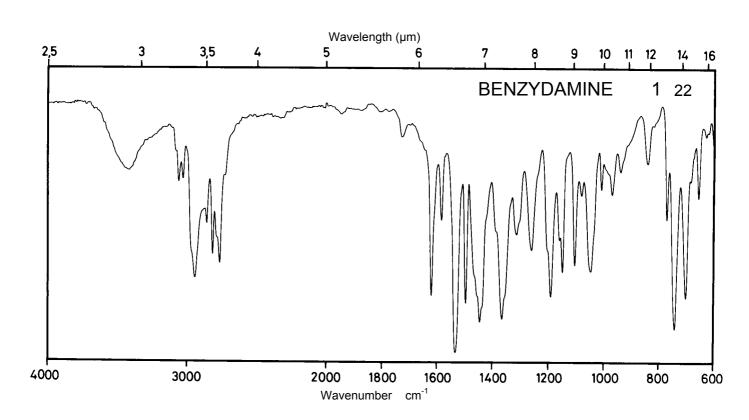





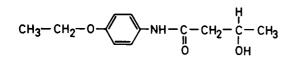


223.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	249 nm		246 nm	246 nm
E ^{1%} 1cm	709		537	537
ε	15830		11990	11990

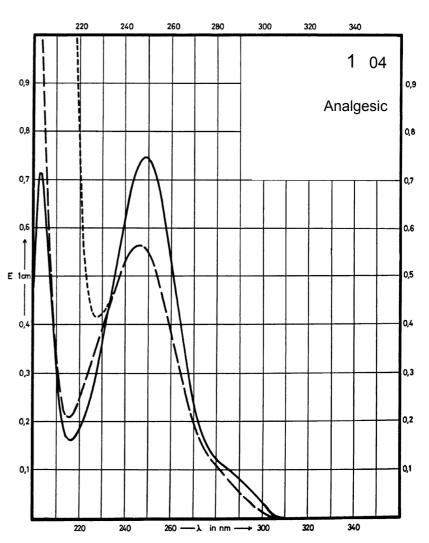

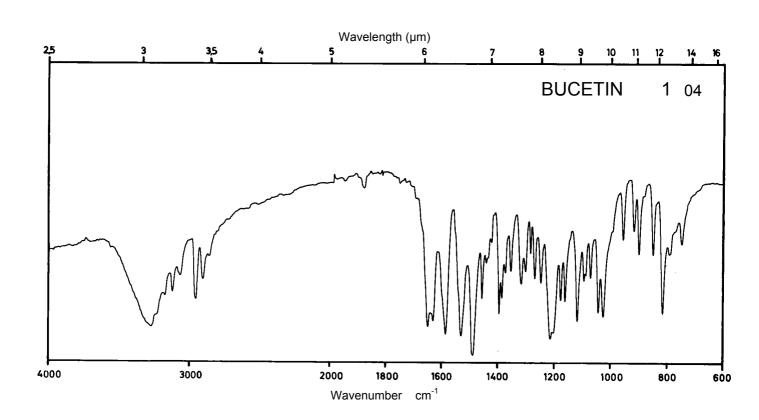





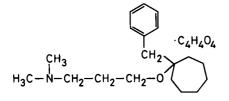

405.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			258 nm	
E ^{1%} 1cm			14.5	
ε			590	

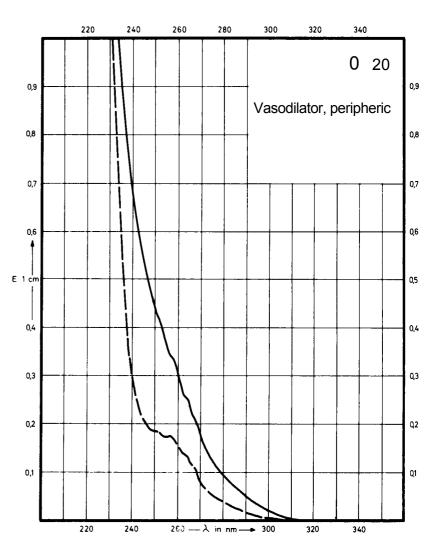

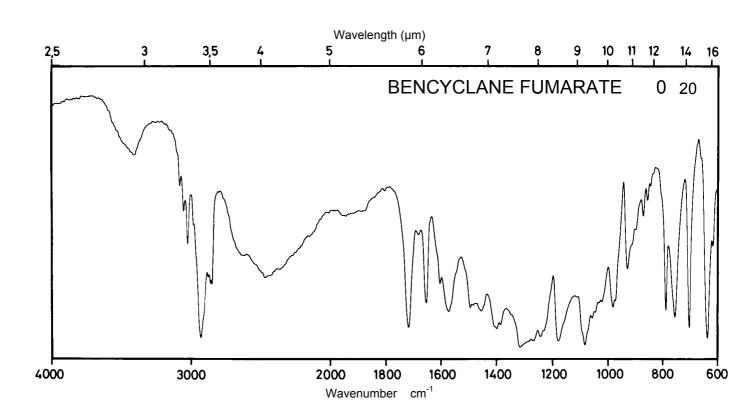







155.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				228 nm
E ^{1%} 1cm				1079
ε				16750

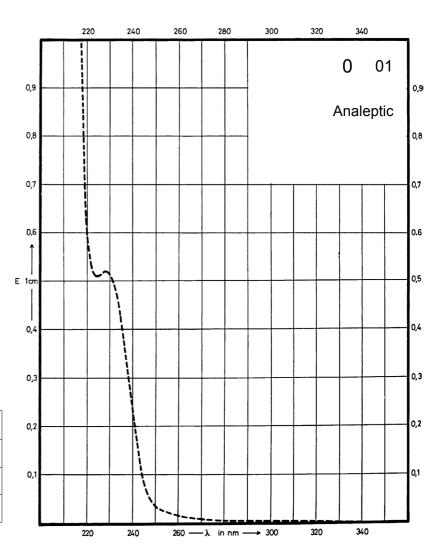

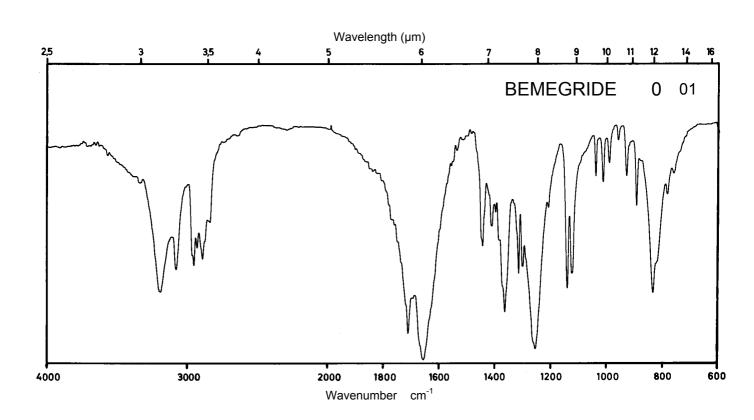


340 260 280 300 320 220 240 **BAMBERMYCINS** Name 4 94 0,9 0,9 Antibiotic 0,8 0,8 0,7 0,7 0,6 0,6 0,5 E 1 cm 0,4 0,4 Concentration 5 mg / 100 ml 0,3 0,3 Solvent Methanol Water 0.1 M HCI 0.1 M NaOH 0,2 Symbol Maximum of 257 nm 245 nm 257 nm absorption 0,1 0,1 64 38 62

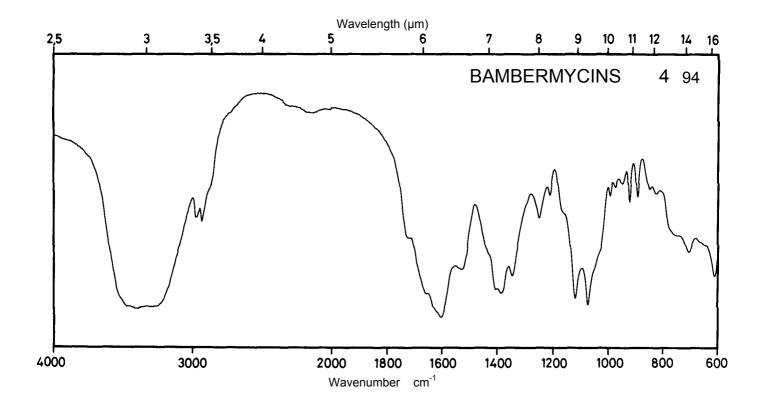

ε



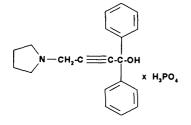



 $\,M_r\,$

389.4


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 258 nm 253 nm			
E ^{1%} 1cm	7.9 10.5 8.9			
3	308 409 347			

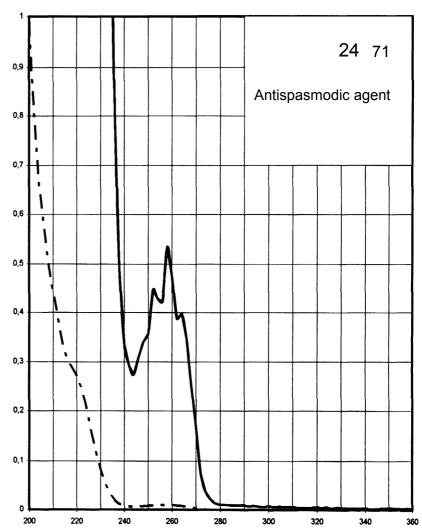

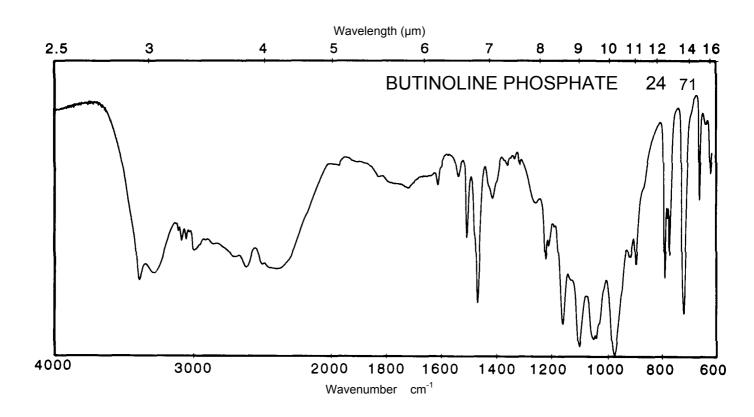





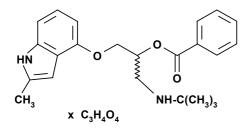


 M_r 484.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm 222 nm	267 nm 220 nm	267 nm 220 nm	267 nm
E ^{1%} 1cm	229 990	199 882	200 885	183
3	11100 48000	9650 42800	9680 42900	8870

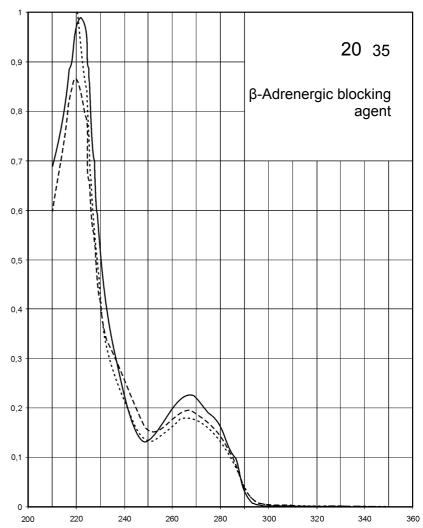

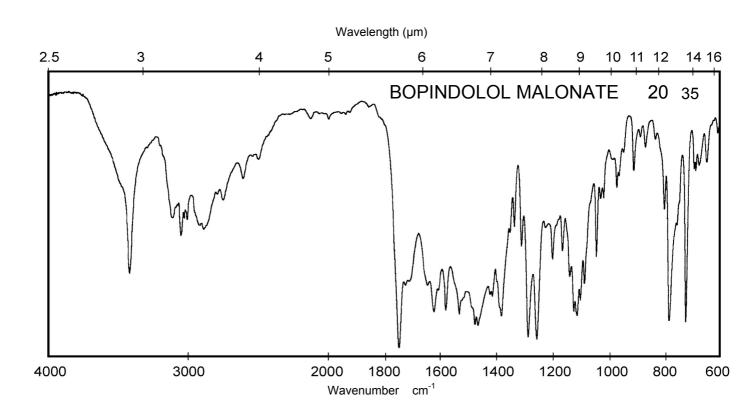




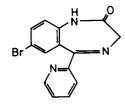


316.2

Concentration 0.5 mg / 100 ml 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	320 nm 233 nm		Decom- position observed	348 nm 237 nm
E ^{1%} 1cm	61 1050			80 920
ε	1930 33200			2530 29090

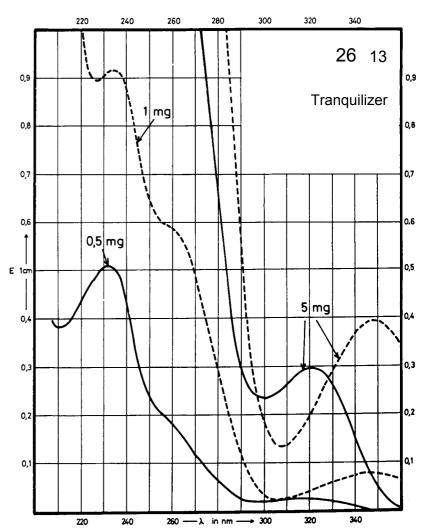


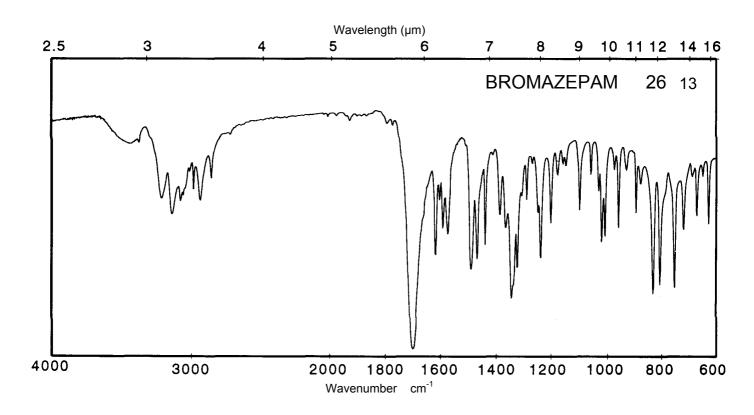





194.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm	273 nm	270 nm	273 nm
E ^{1%} 1cm	475	515	495	510
ε	9220	10000	9610	9900

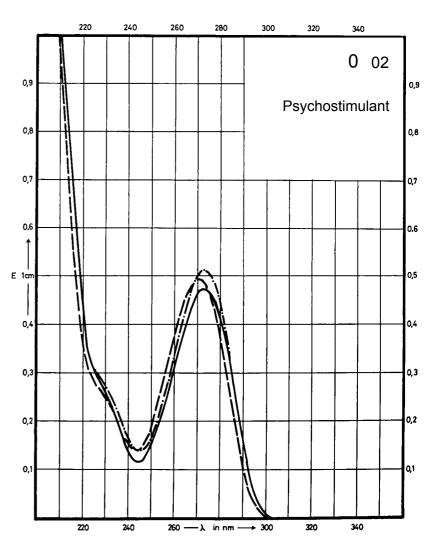

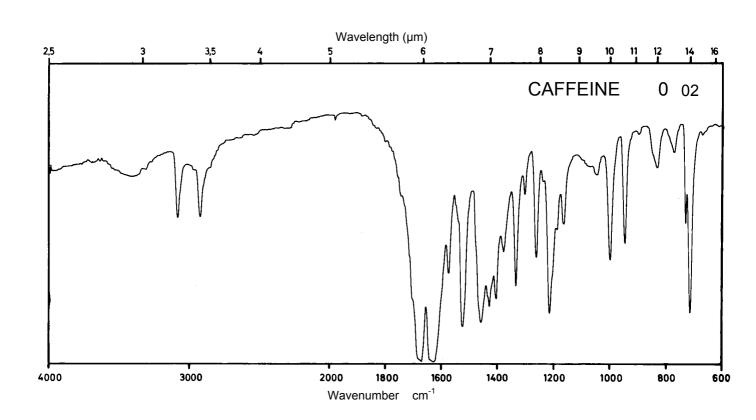





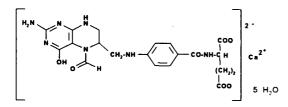



601.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 221 nm	286 nm 220 nm	Decom- position observed	221 nm
E ^{1%} 1cm	560 595	512 571		507
ε	33700 35800	30800 34400		30500

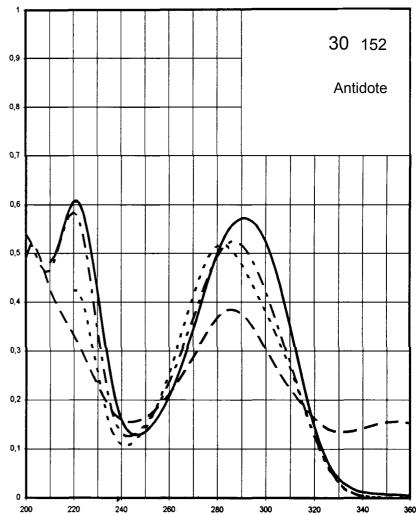

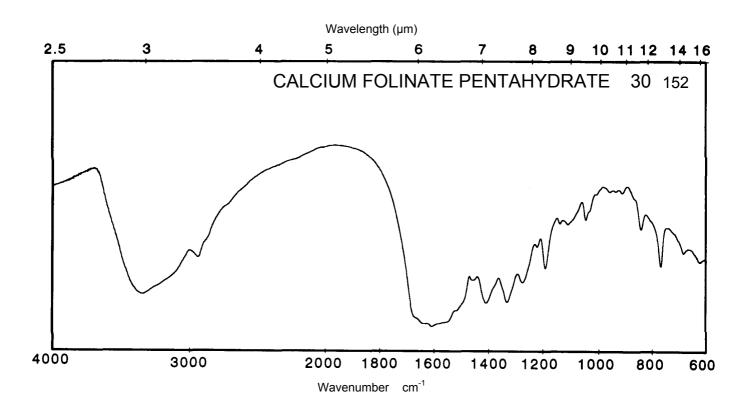


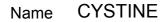


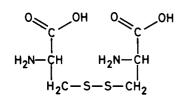



240.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

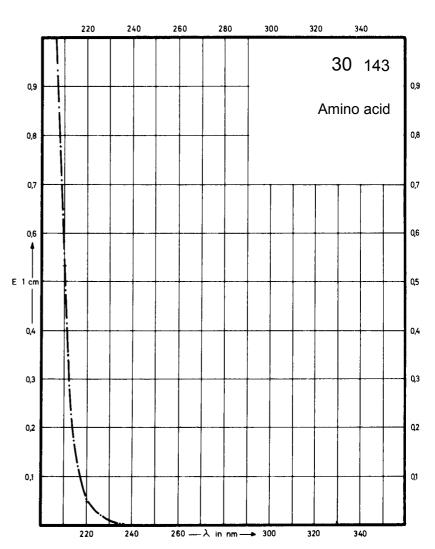


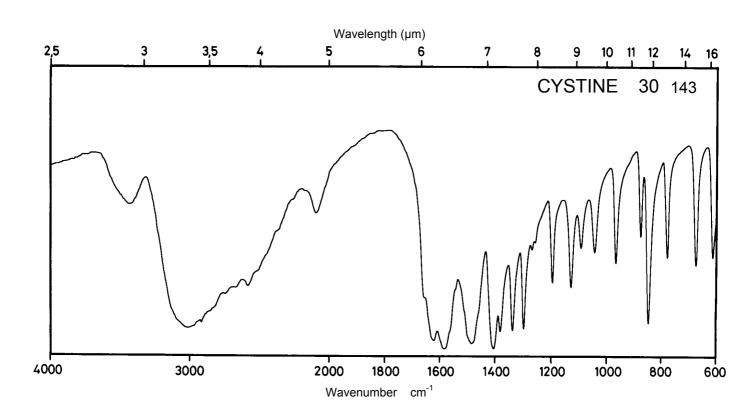

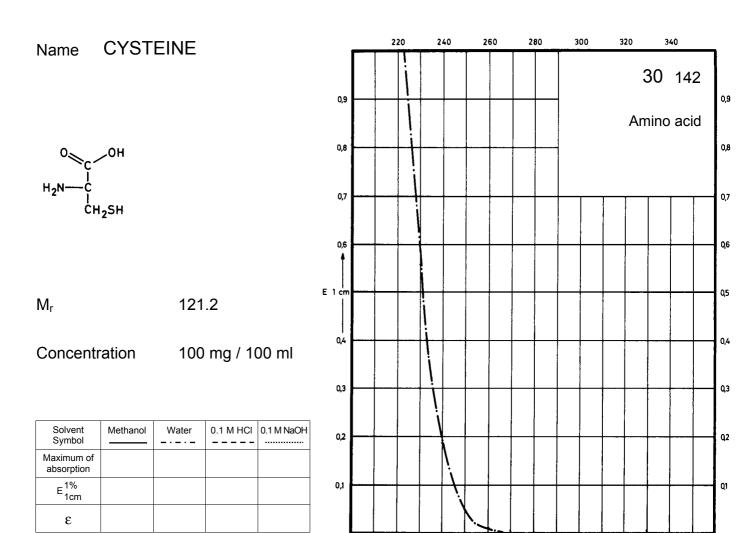

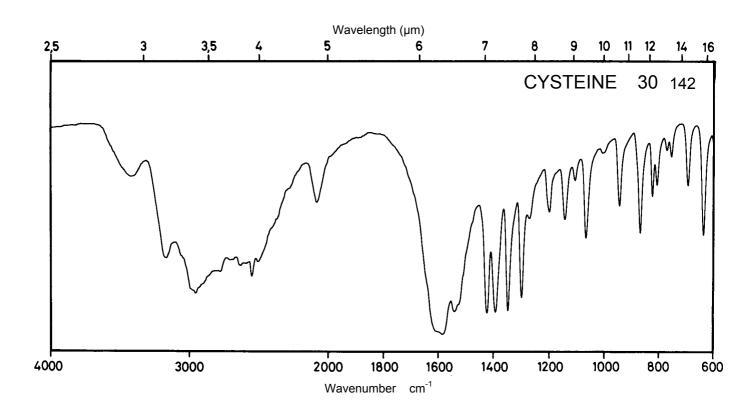





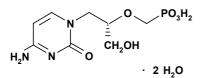


315.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	286 nm	282 nm 213 nm	282 nm 214 nm	274 nm
E ^{1%} 1cm	423	403 311	413 308	283
3	13300	12700 9800	13000 9700	8930

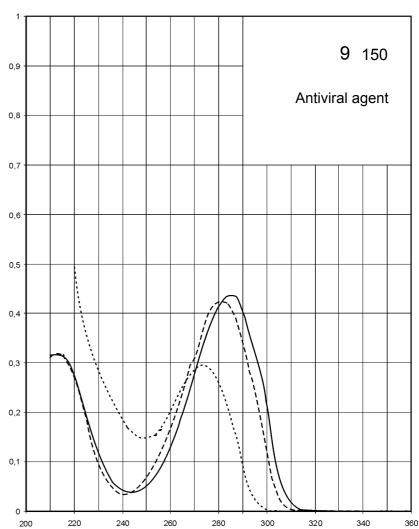

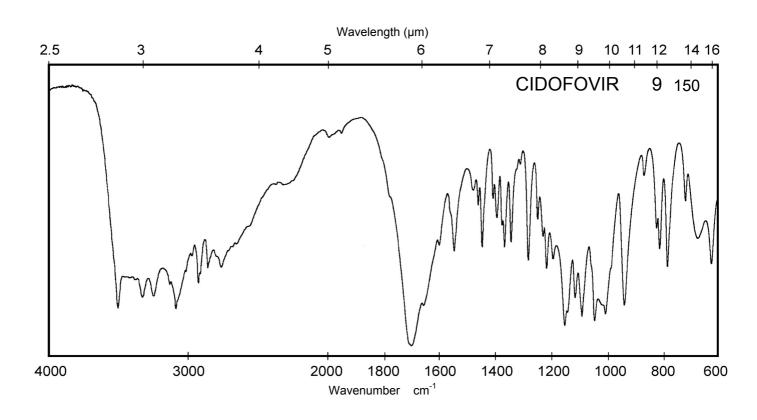


Name CROCOCONAZOLE HYDROCHLORIDE

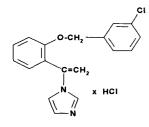



 M_{r}

347.3


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	286 nm	284 nm	285 nm	279 nm
E ^{1%} 1cm	88	76	82	75
ε	3050	2630	2840	2600

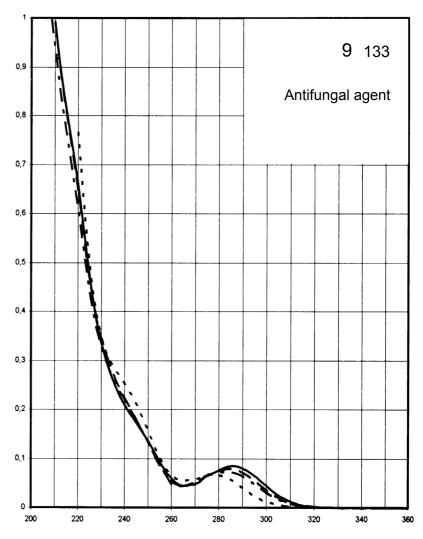


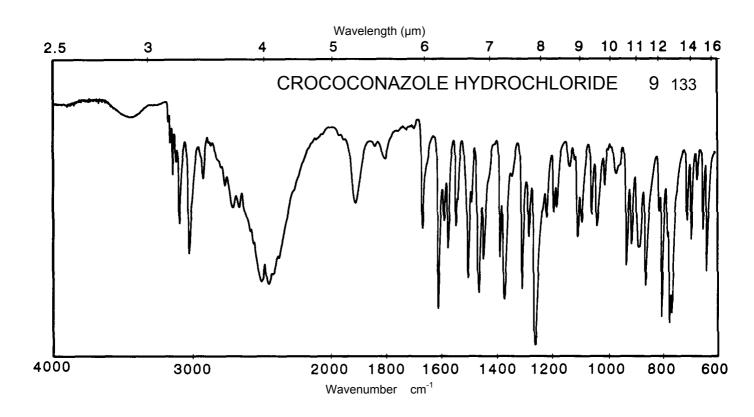

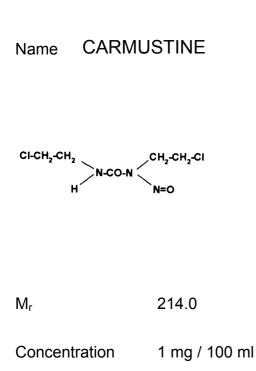


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	231 nm	229 nm	229 nm	Decom- position observed
E ^{1%} 1cm	278	271	269	
ε	6000	5800	5700	

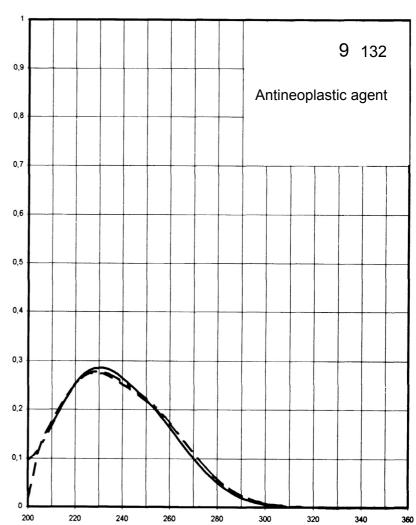

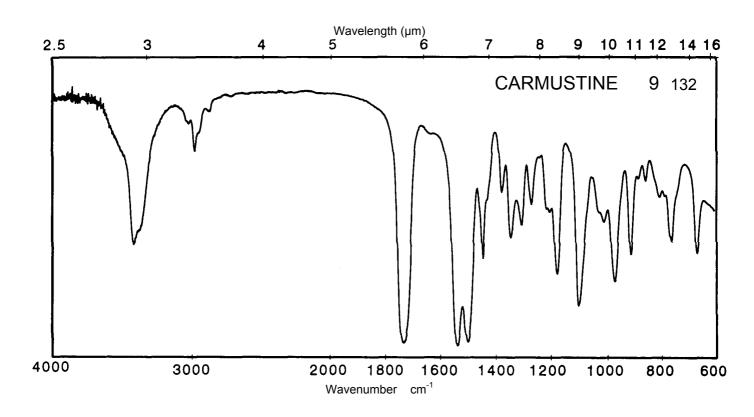



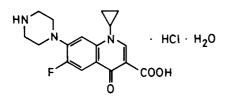



385.8


Concentration


0.55 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	330 nm 317 nm 280 nm		329 nm 315 nm 277 nm	335 nm 323 nm 271 nm
E ^{1%} 1cm	302 316 1036		283 318 1076	378 378 885
ε	11700 12200 40000		11000 12260 41500	14600 14600 34100

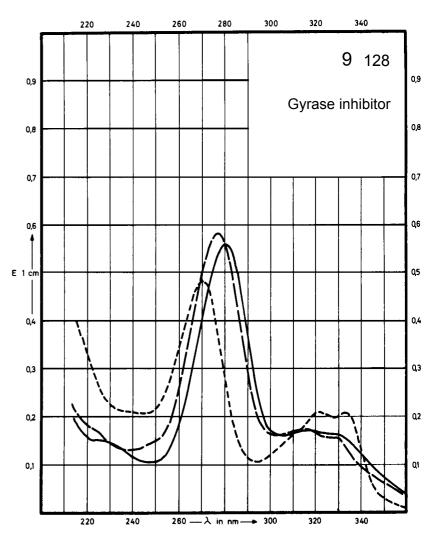

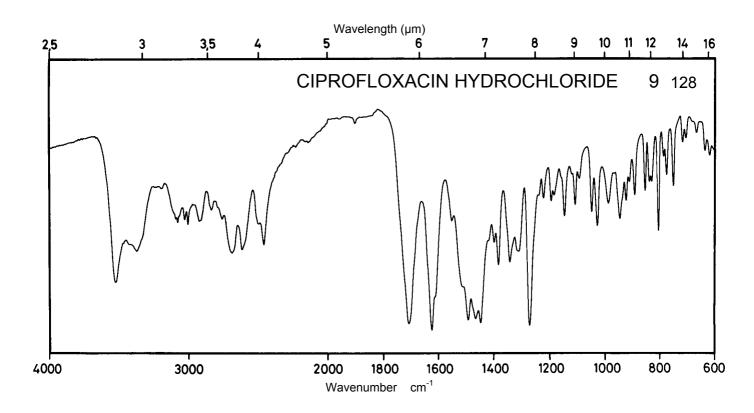

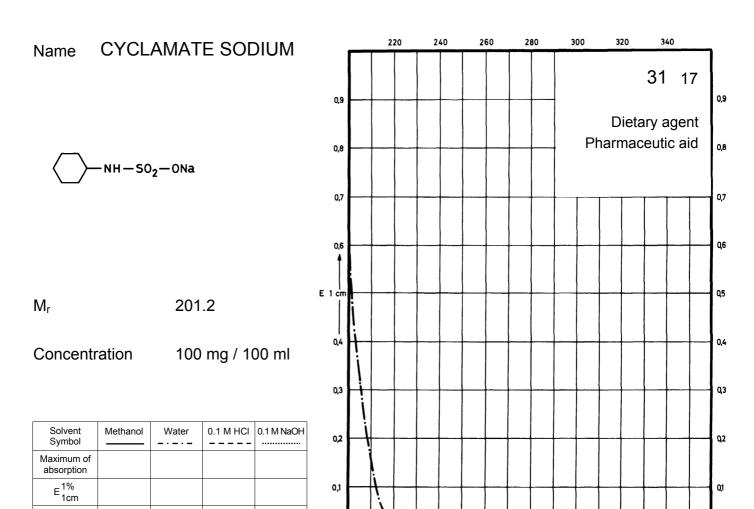

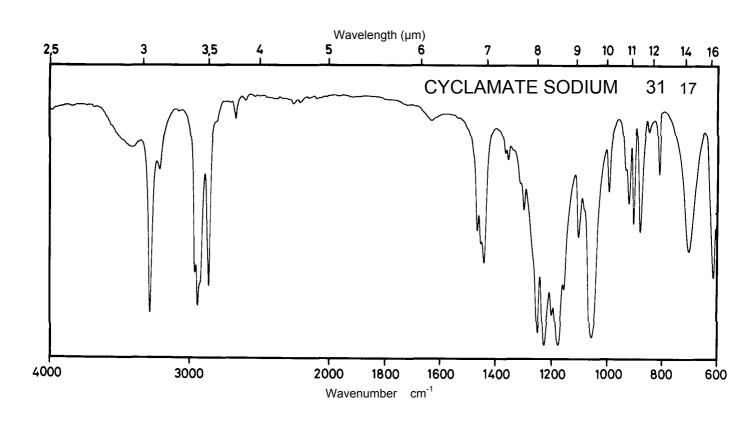


240


260 -

—λ in nm —→ 300

340

3

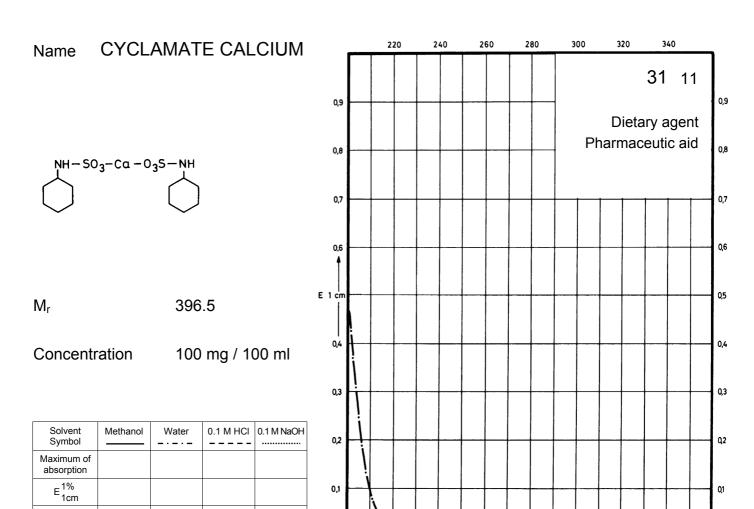

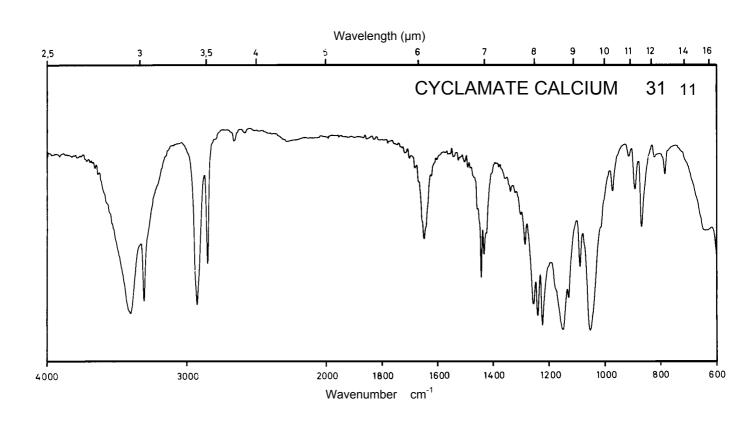


240

260 — λ in nm → 300

340

3

Name DEMBREXINE HYDROCHLORIDE

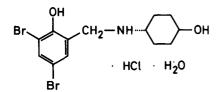



 M_r

433.6


Concentration 2 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	294 nm		291 nm	310 nm 249 nm
E ^{1%} 1cm	72		65	111 243
ε	3130		2800	4800 10500

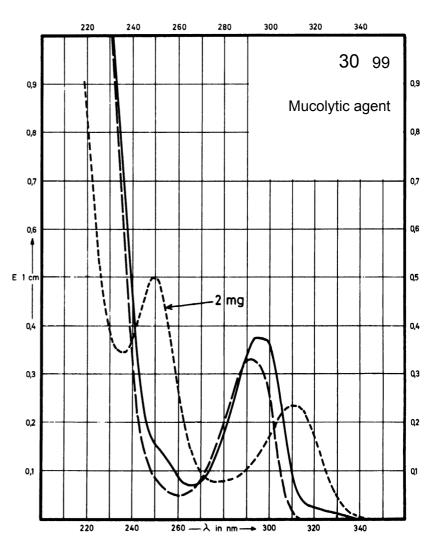

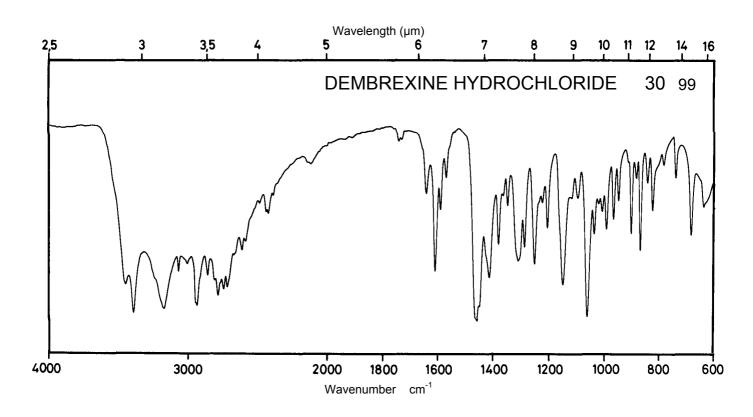





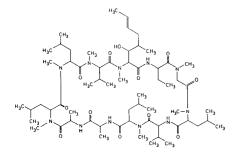


 M_r 1202.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

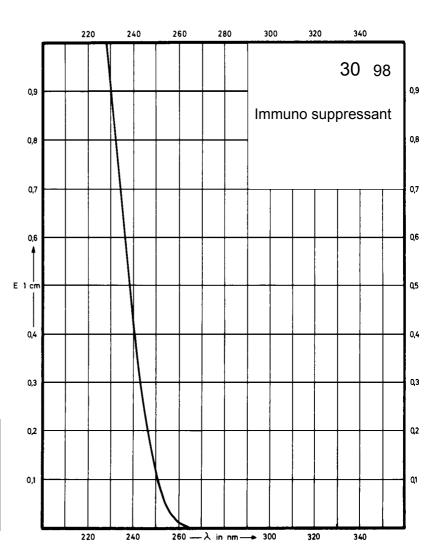

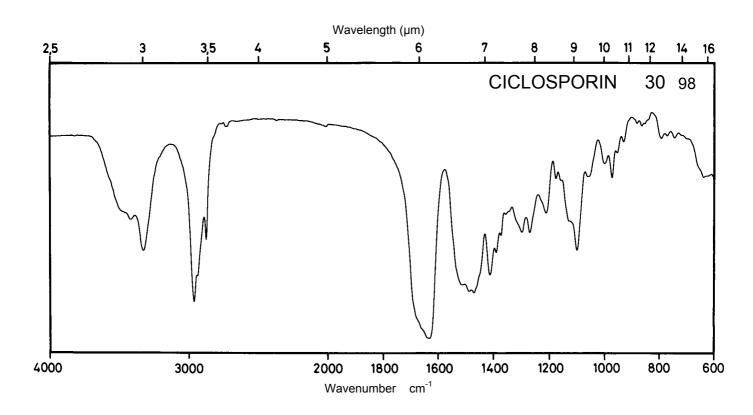





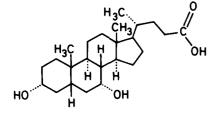


 M_{r} 392.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

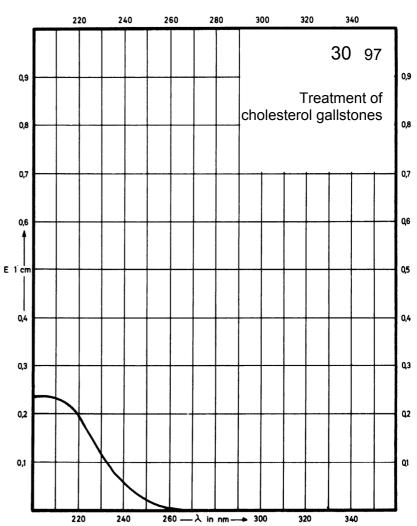

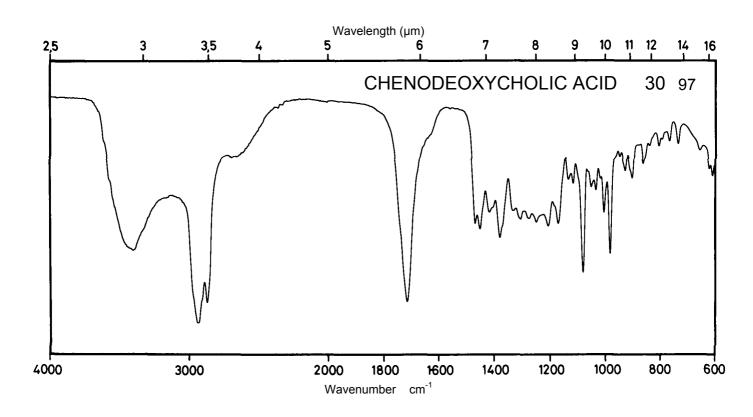


Name CYTIDINE DISODIUM PHOSPHATE

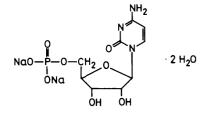


 M_{r}


403.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		280 nm	271 nm
E ^{1%} 1cm	186		288	200
ε	7500		11600	8000

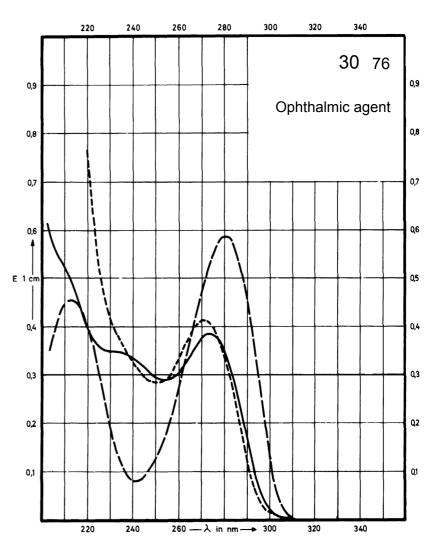

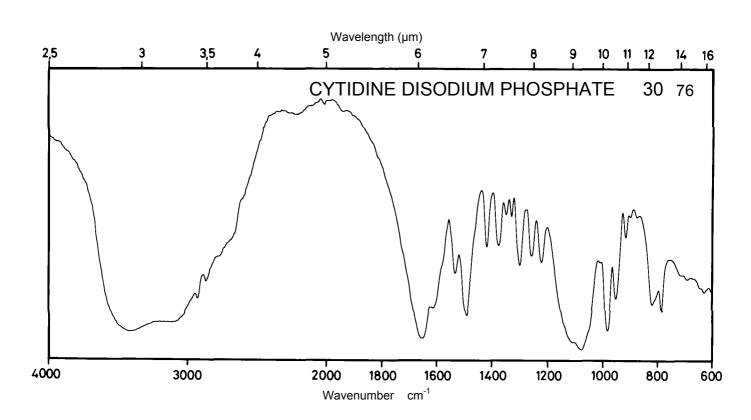





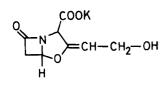


237.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm			
E ^{1%} 1cm	Decom- position observed	Decom- position observed	Decom- position observed	Decom- position observed
3				

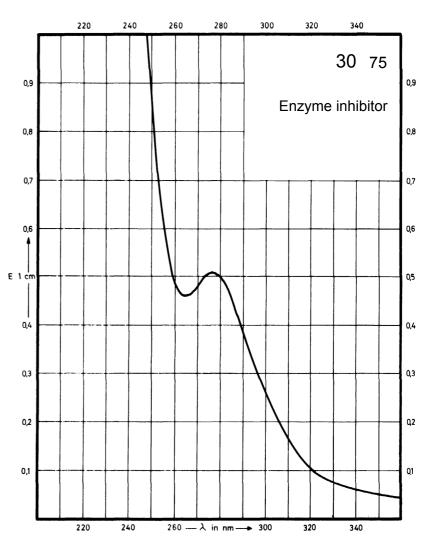

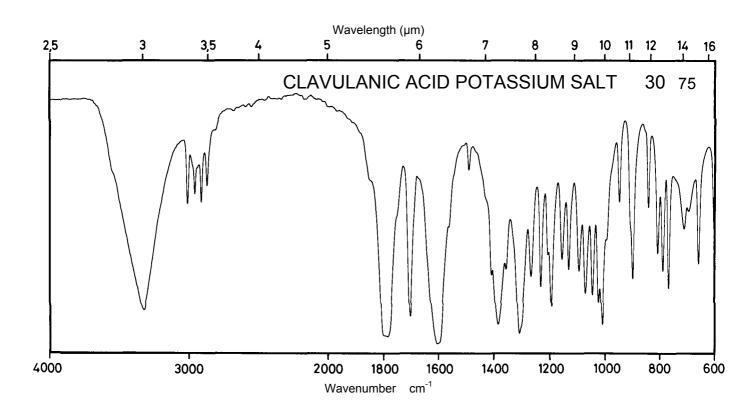




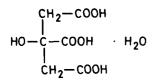


 $\,M_r\,$


210.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

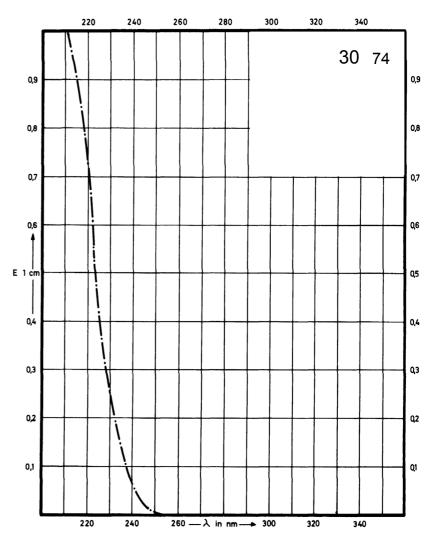

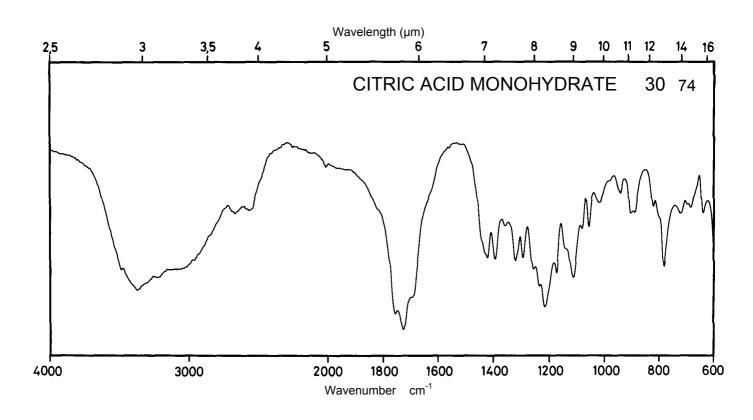



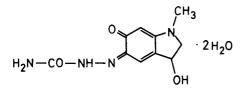



272.3

Concentration


0.75 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	440 nm 356 nm		376 nm	442 nm
E ^{1%} _{1cm}	380 970		1170	1320
ε	10400 26400		31900	35900

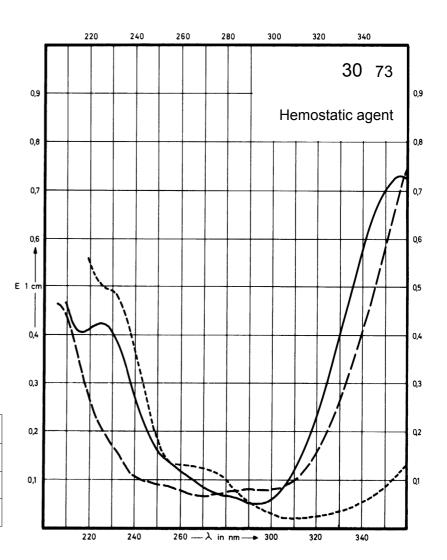

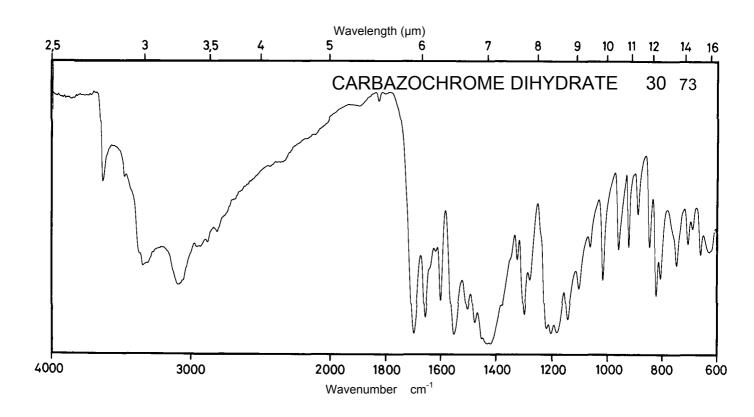




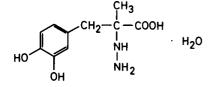


 M_r 244.2


Concentration 2 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm	281 nm	281 nm	Decom- position observed
E ^{1%} 1cm	129	114	114	
ε	3150	2790	2790	

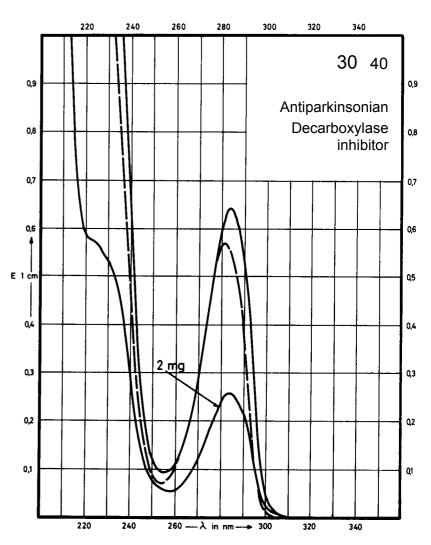

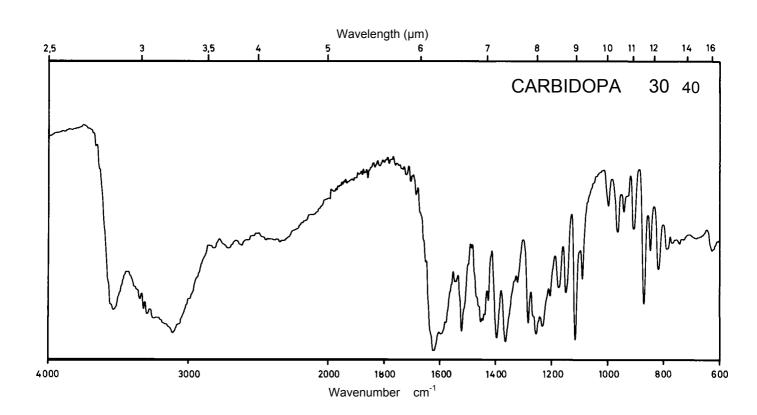




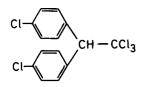


354.5

Concentration 1.25 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm 236 nm			
E ^{1%} 1cm	19 510			
ε	674 18080			

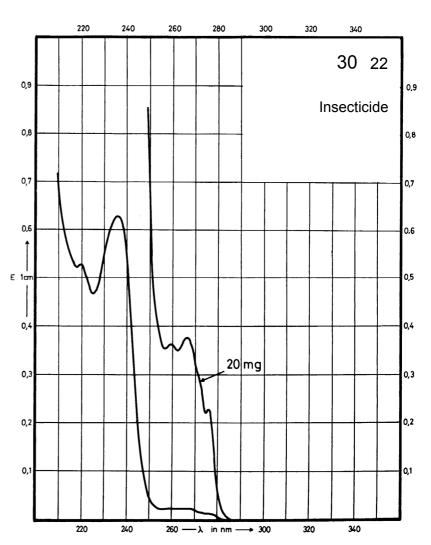

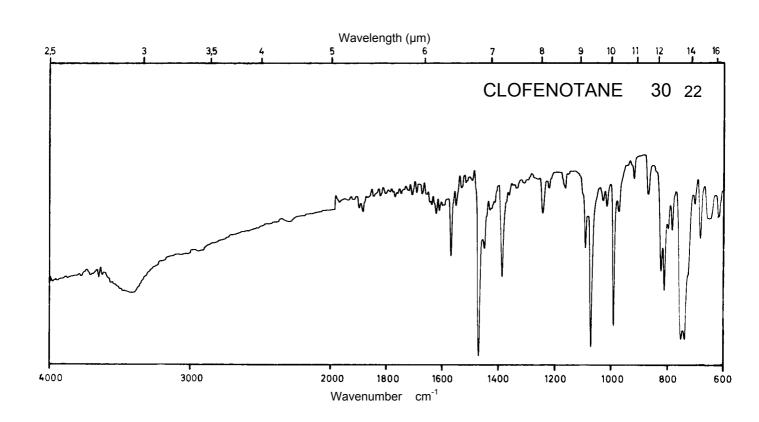



Name 4-CHLORO-5-SULFAMOYL ANTHRANILIC ACID

 M_{r}


250.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	329 nm 268 nm 229 nm		330 nm 268 nm 230 nm	318 nm 262 nm
E ^{1%} 1cm	170 675 1740		174 622 1940	130 540
ε	4260 16930 43650		4370 15600 48730	3260 13520

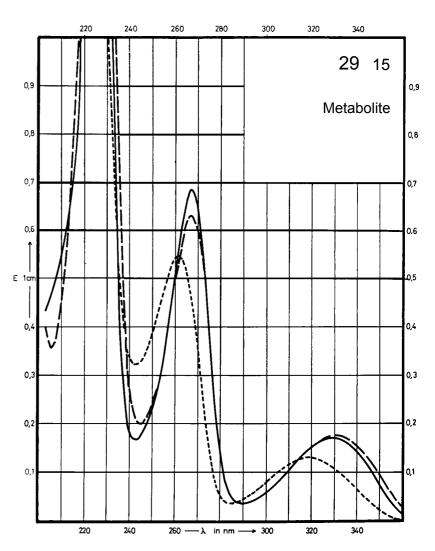

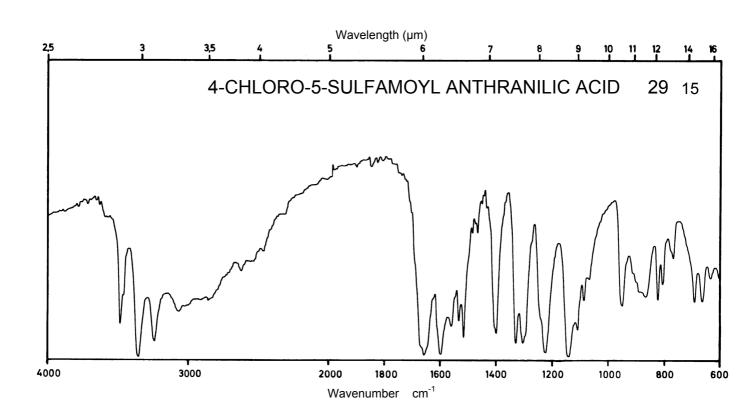





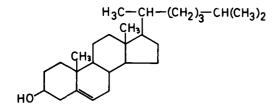


386.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

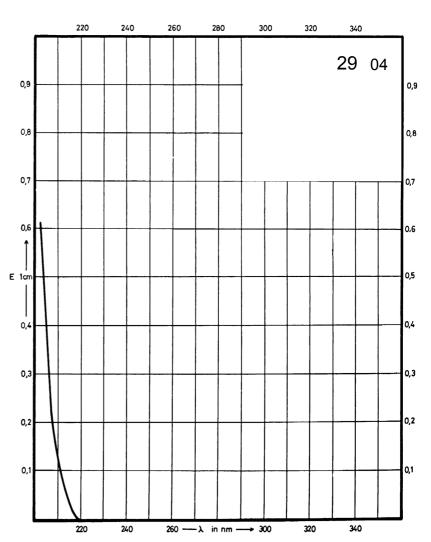

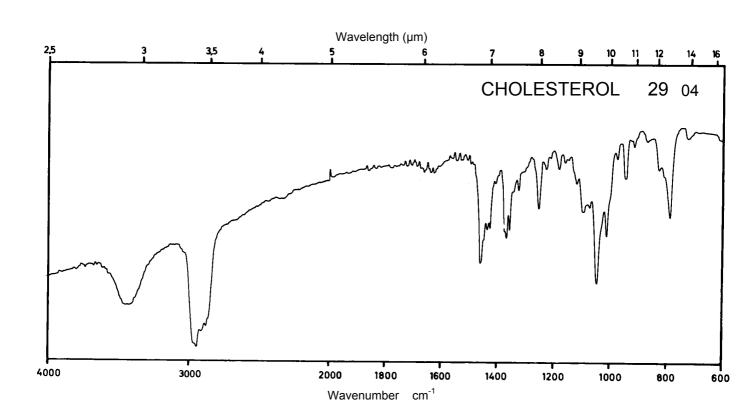





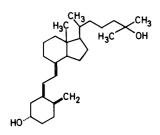


 M_r 400.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm			
E ^{1%} 1cm	440			
3	17600			

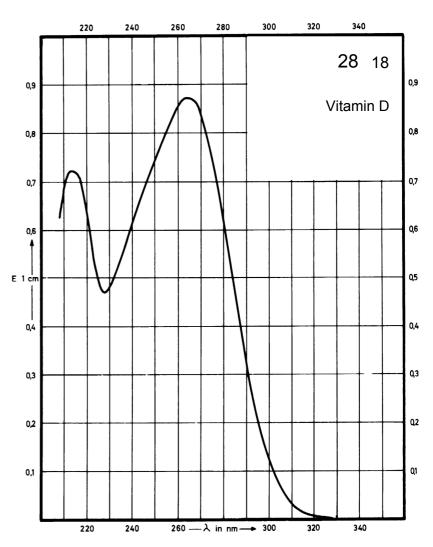

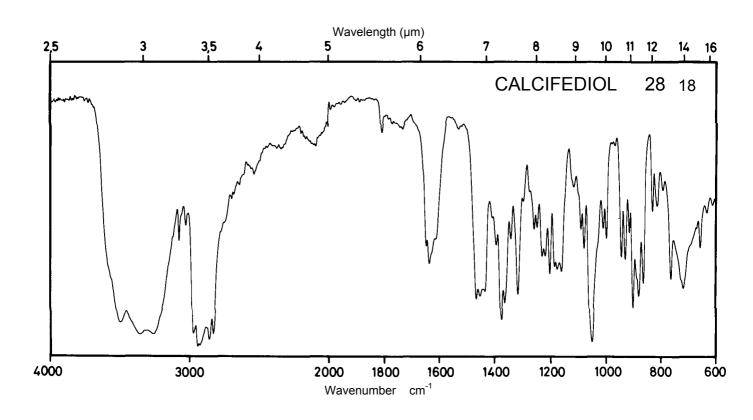




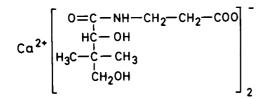


 M_r


476.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

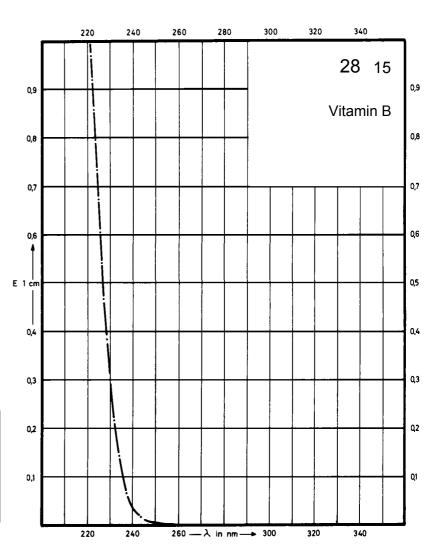

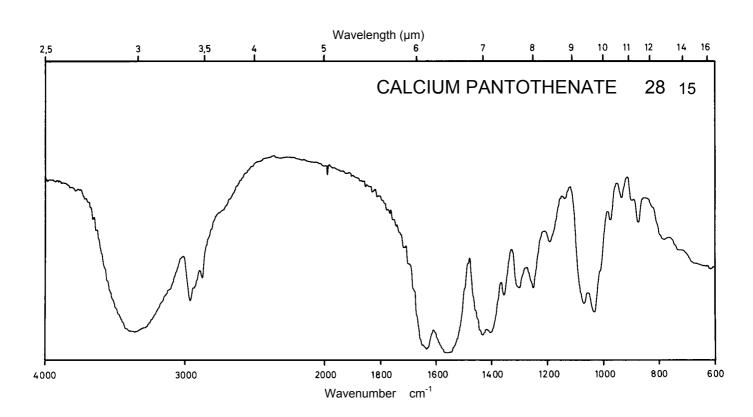





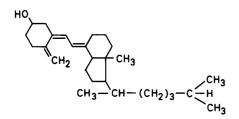

384.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm			
E ^{1%} 1cm	479			
3	18430			

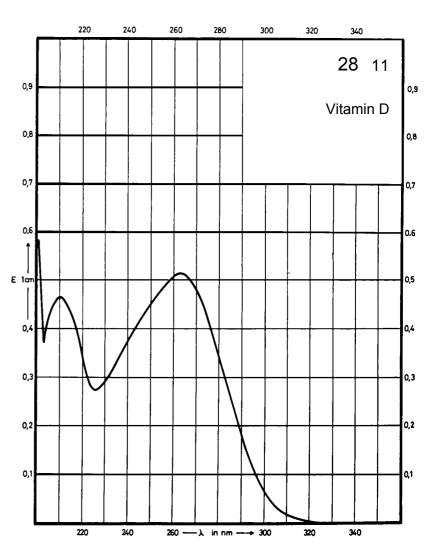

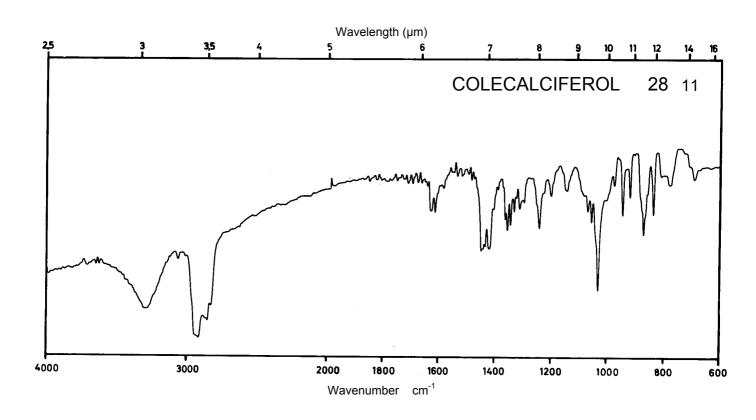



220 240 260 280 300 320 340

28 02

0,9

0,8

0,7


0,6


1 tem

 M_r


1355.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		360 nm 277 nm		
E ^{1%} 1cm		188 110		
3		25480 14910		

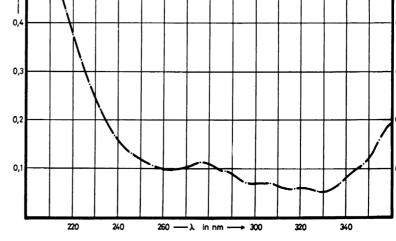

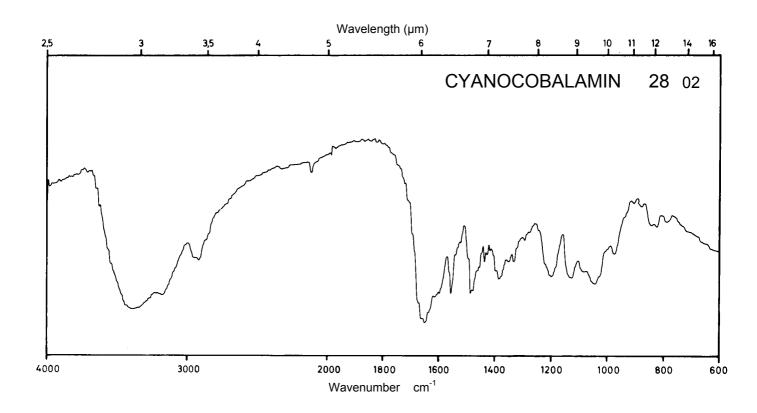





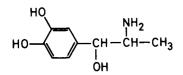


183.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			279 nm	
E ^{1%} 1cm			146	
ε			2670	

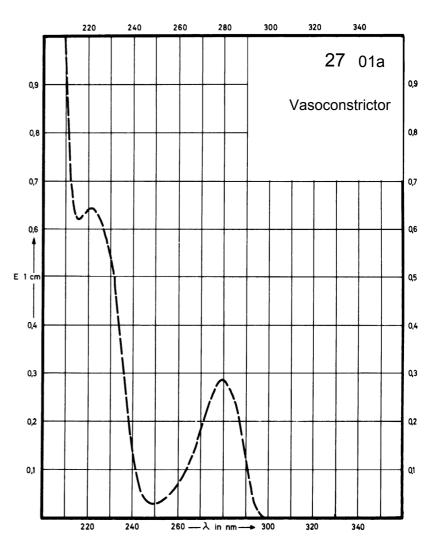

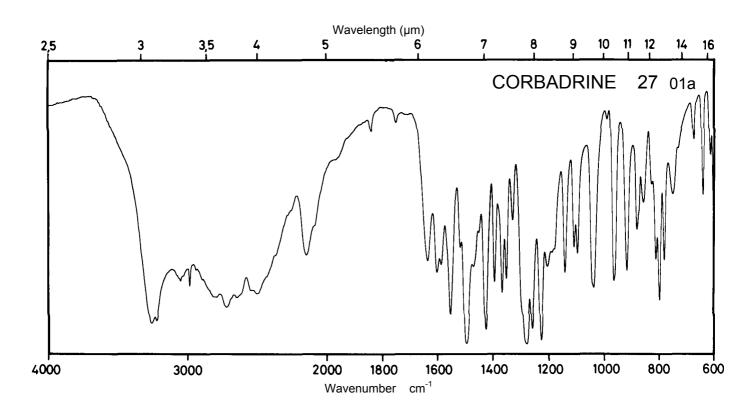


Name CORBADRINE HYDROCHLORIDE

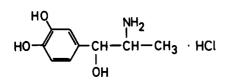


 M_{r}


219.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm		279 nm	
E ^{1%} 1cm	142		122	
3	3120		2680	

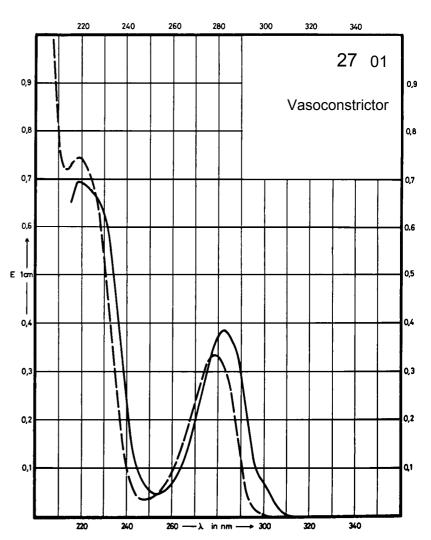

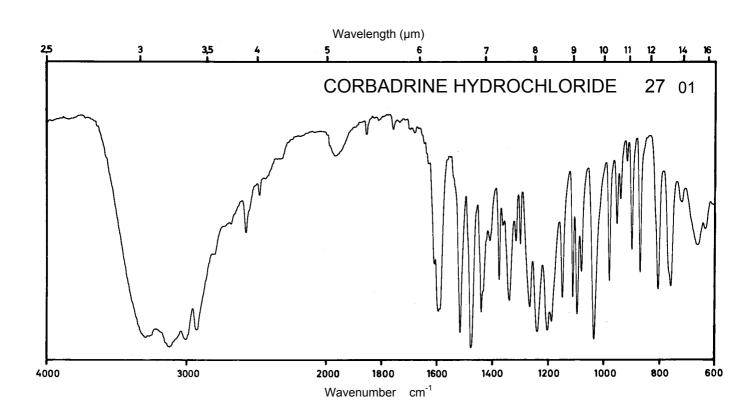





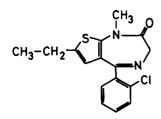


318.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	320 nm 243 nm	319 nm 244 nm	393 nm 261 nm	Decom- position observed
E ^{1%} 1cm	87 680	86 680	148 623	
3	2770 21680	2740 21680	4720 19860	

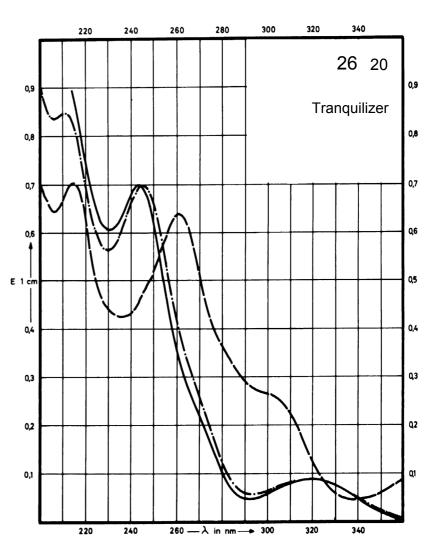

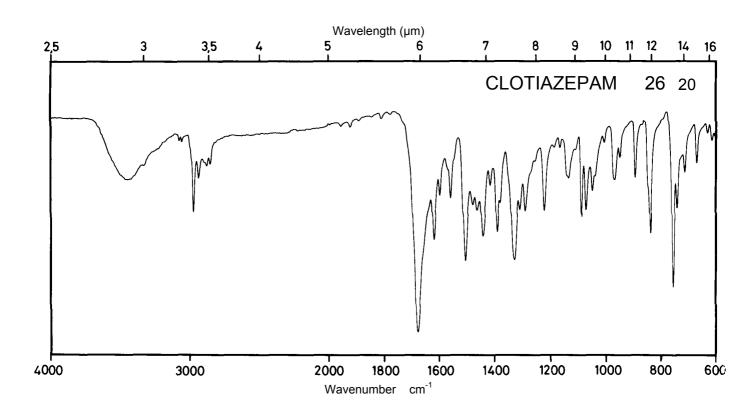




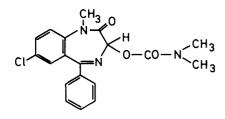


 M_r 371.8


Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm 231 nm		Decom- position observed	Decom- position observed
E ^{1%} 1cm	67 937			
ε	2490 34840			

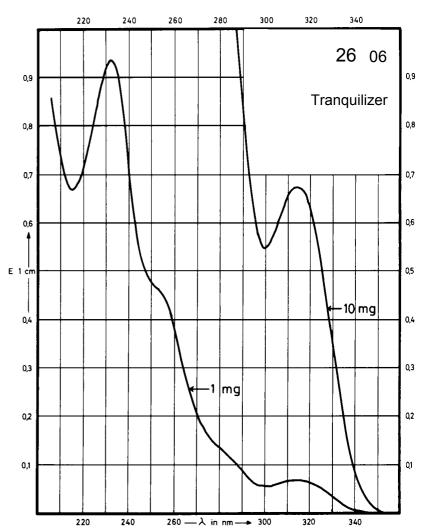

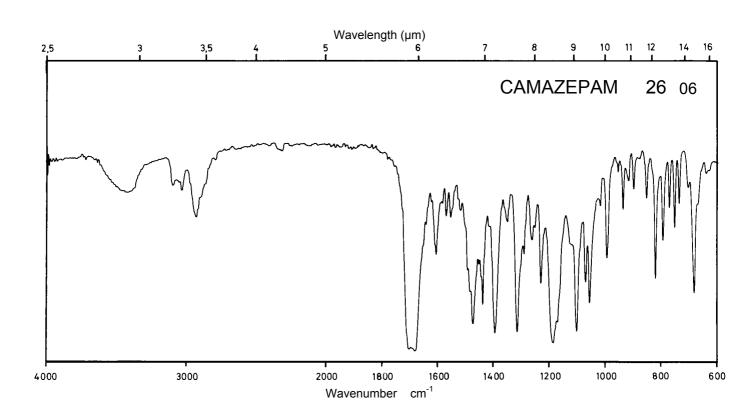




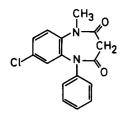


300.7


Concentration 0.5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	290 nm 230 nm		289 nm 230 nm	286 nm
E 1% 1cm	76 1385		76 1373	193
ε	2290 41650		2290 41290	5800

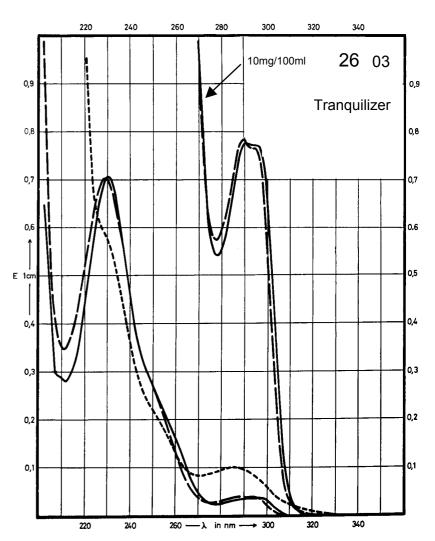

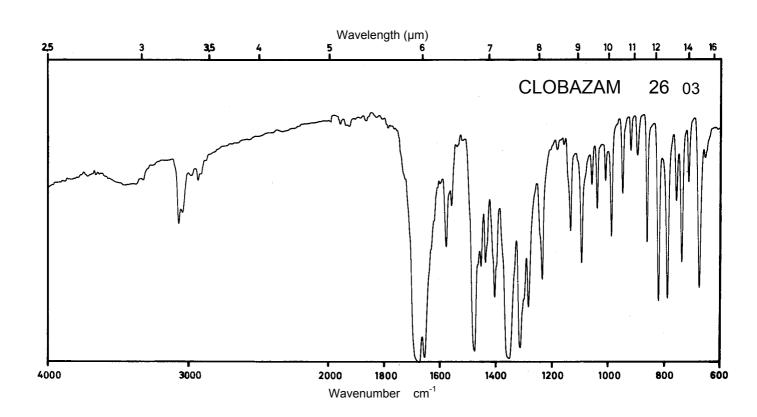

Name CHLORDIAZEPOXIDE HYDROCHLORIDE

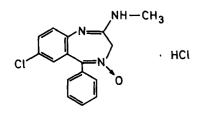


 M_{r}


336.2

Concentration


0.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm 245 nm		307 nm 245 nm	260 nm
E ^{1%} 1cm	255 945		290 1025	1025
ε	8600 31800		9750 34500	34500

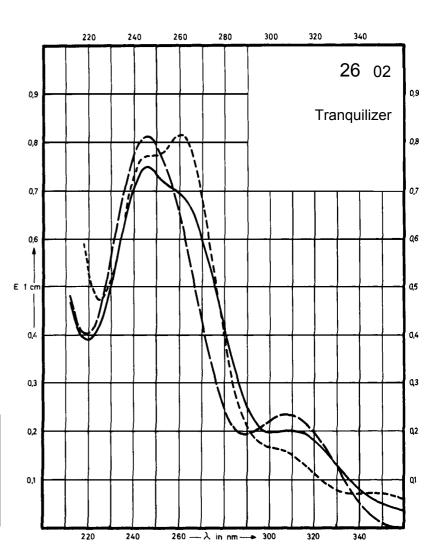

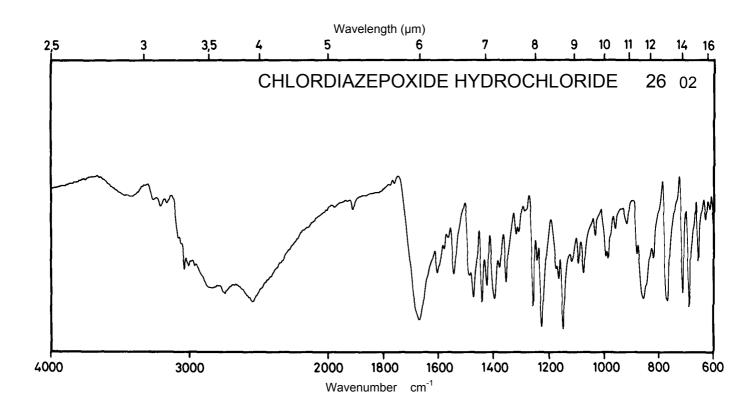


Name CLORAZEPATE, DIPOTASSIUM

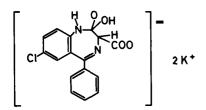



 M_r 408.9

Concentration 1 mg / 100 ml 11 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 225 nm		360 nm 283 nm 237 nm	344 nm
E ^{1%} 1cm	53 824		107 335 747	58
ε	2170 33690		4380 13700 30540	2370

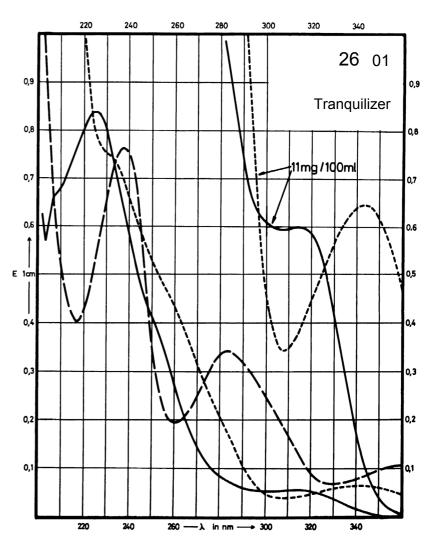

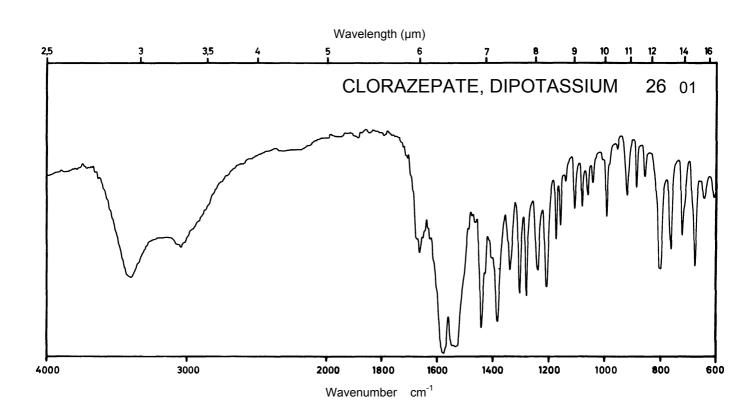


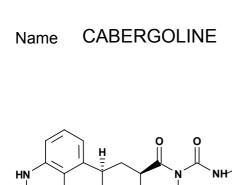


451.7

N(CH₃)₂

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm 223 nm	280 nm 222 nm	279 nm 221 nm	280 nm
E ^{1%} 1cm	162 879	144 814	151 834	141
3	7320 39700	6520 36800	6840 37700	6380

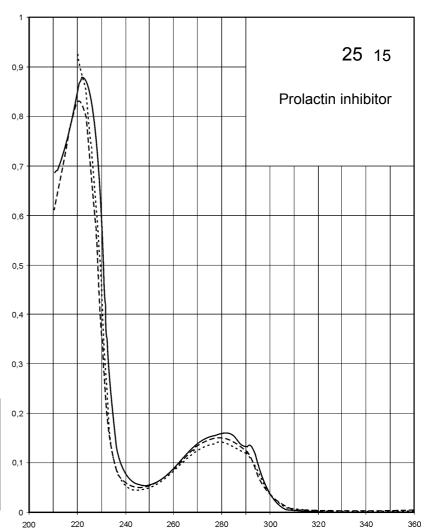

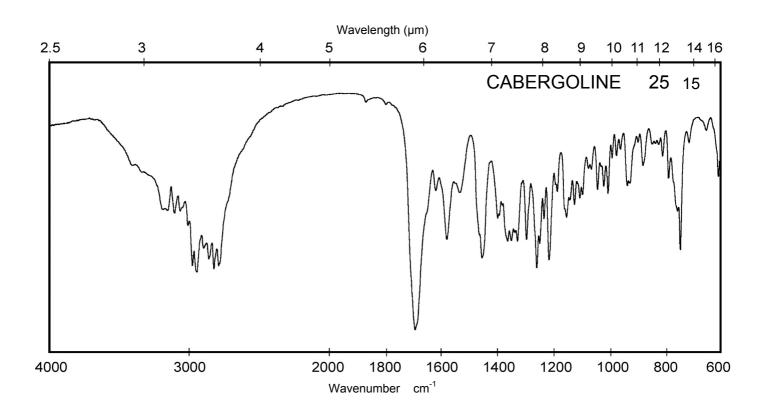





H - Cys - Gly - Asn - Leu - Ser - Thr - Cys - Met - Leu - Gly - Thr - Tyr - Thr - Gln - Asp - Phe - Asn - Lys - Phe - His - Thr - Phe - Pro - Gln - Thr - Ala - Ile - Gly - Yal - Gly - Ala - Pro - NH₂


 M_r 3417.9

	olvent ymbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
_	imum of orption	278 nm		275 nm	
E	1% 1cm	4.9		4.4	
	ε	1670		1500	

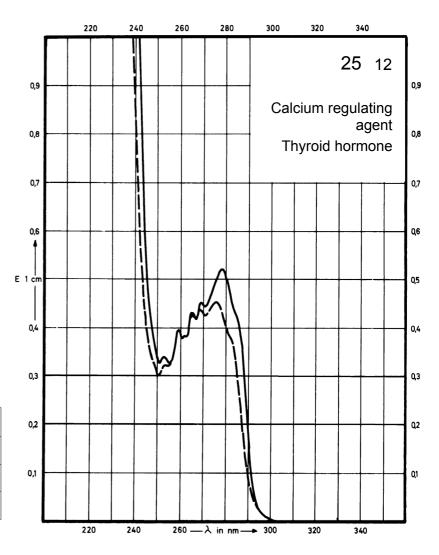

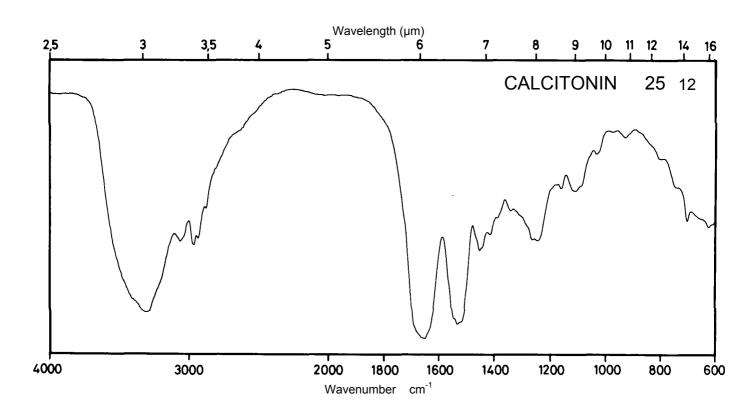





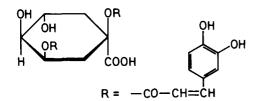


516.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	326 nm 246 nm		321 nm	Decom- position observed
E ^{1%} 1cm	640 435		605	
ε	33000 22500		31200	

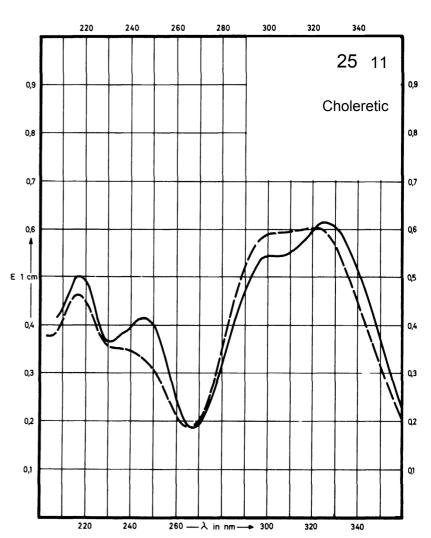

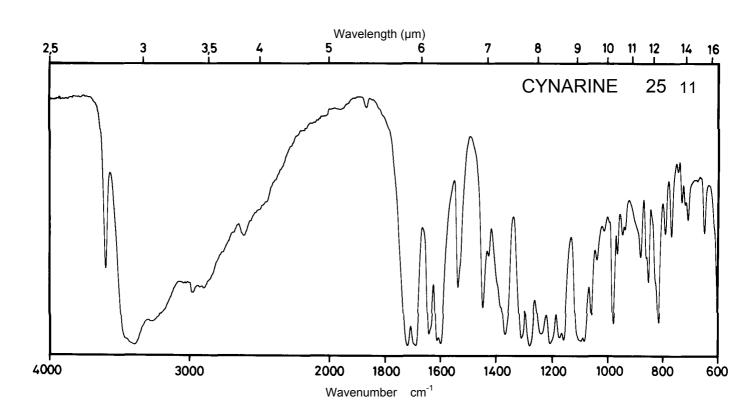





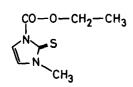


 M_{r} 186.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	298 nm			
E ^{1%} 1cm	Decom- position observed			
ε				

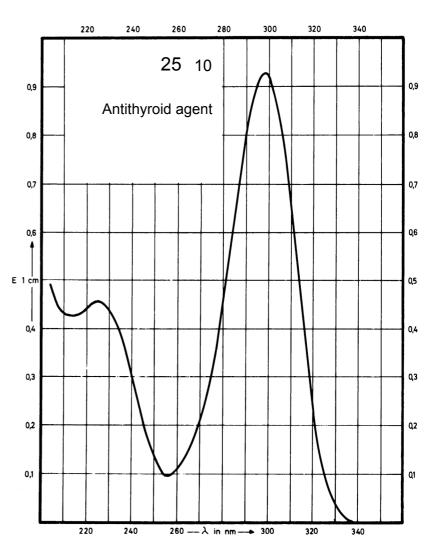

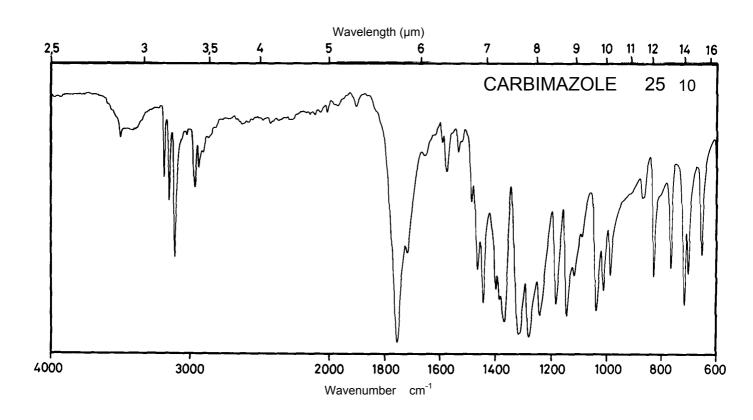





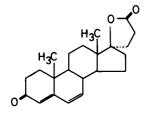


340.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm		292 nm	293 nm
E ^{1%} 1cm	798		770	768
ε	27130		26210	26140

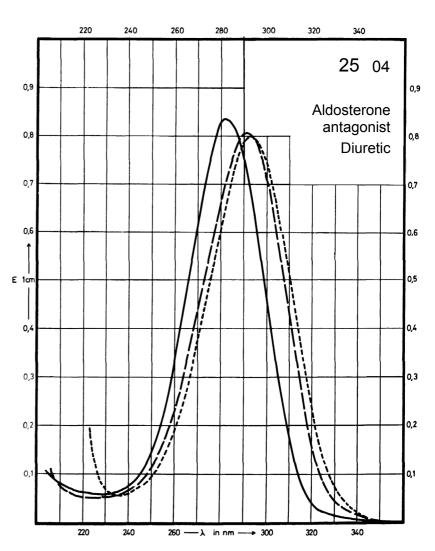

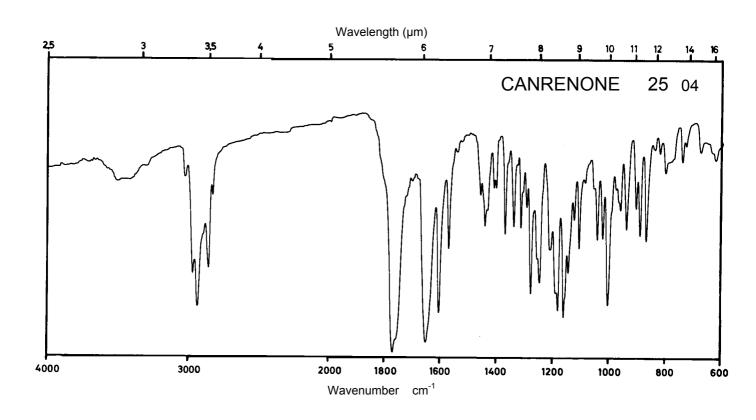





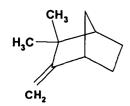



136.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

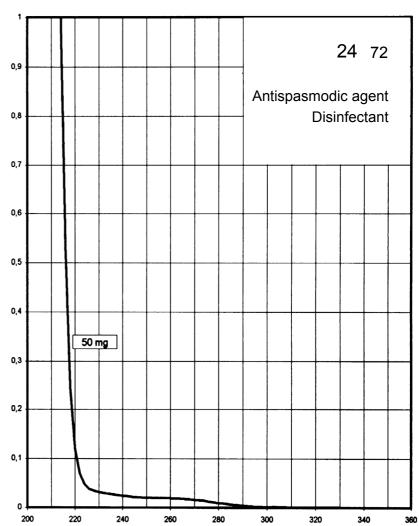


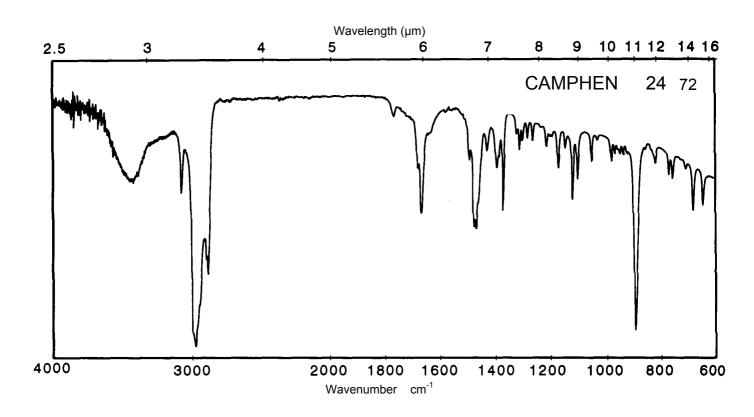

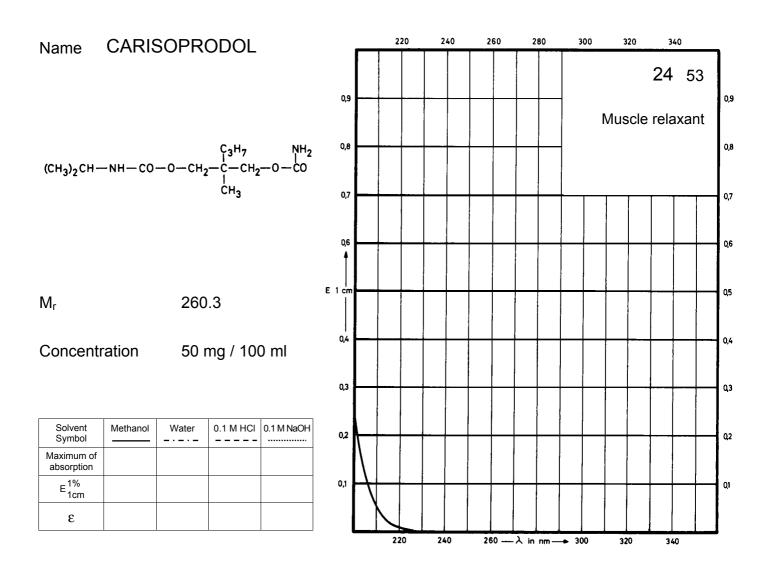

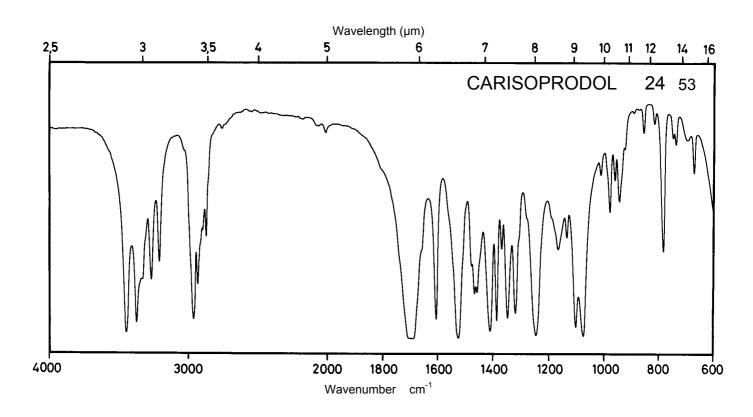





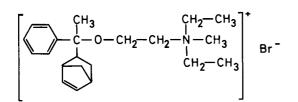


 M_r 408.4

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	264 nm 258 nm		Decom- position	
absorption	252 nm		observed	
E ^{1%} 1cm	3.86 4.93			
TOTT	3.80 158			
3	200			

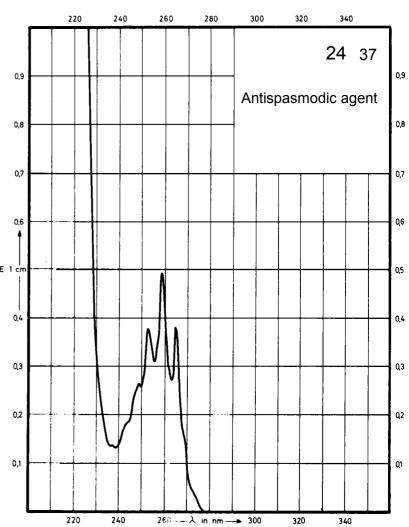

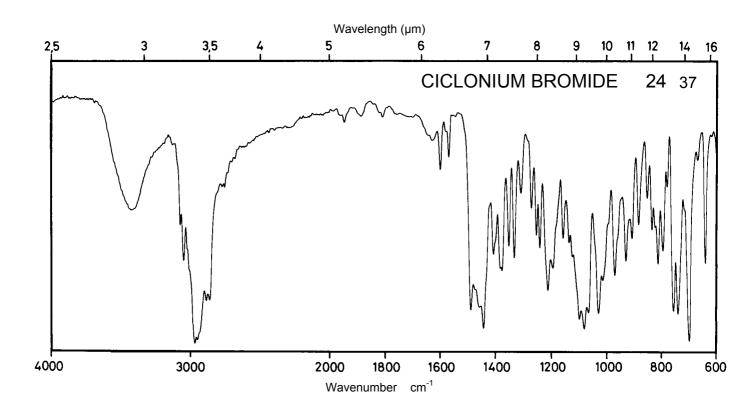





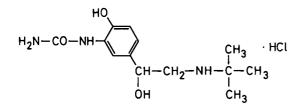


 M_r 303.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm 240 nm		282 nm	300 nm 245 nm
E ^{1%} 1cm	145 390		95	169 359
ε	4400 11800		2900	5150 10900

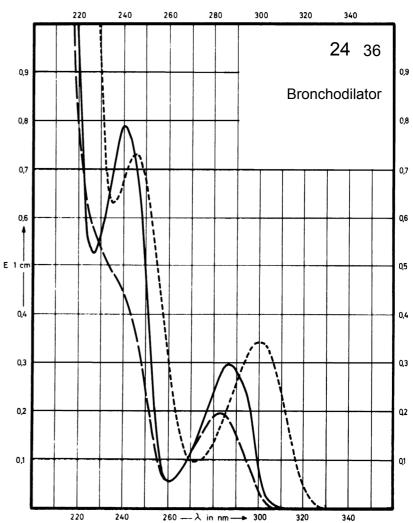

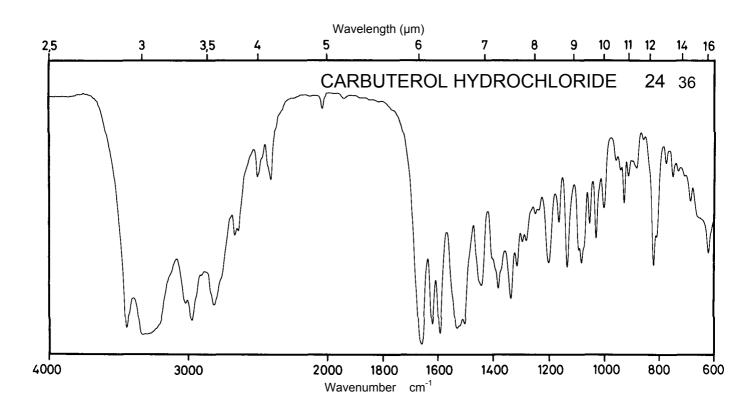




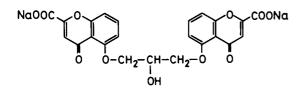


512.3

 M_r

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	320 nm 229 nm	327 nm 238 nm		327 nm 229 nm
E ^{1%} 1cm	210 702	145 540		147 534
3	10760 35960	7430 27660		7530 27360

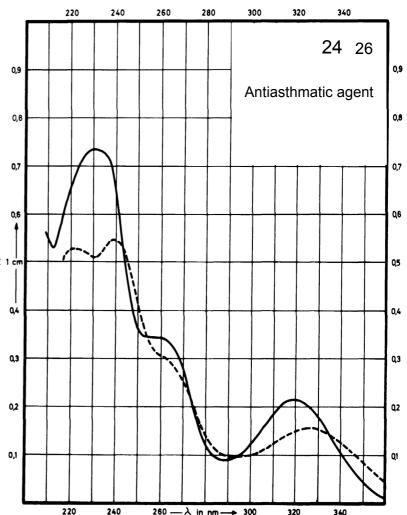

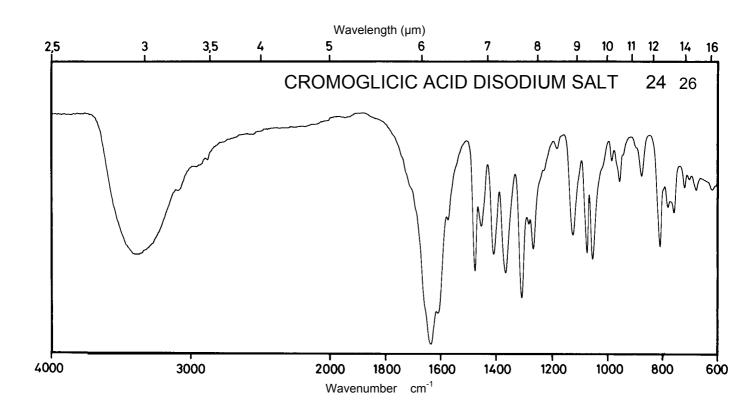

Name CLENBUTEROL HYDROCHLORIDE

$$\begin{array}{c|c} & \text{CH}_3 \\ \text{CH}-\text{CH}_2-\text{NH}-\text{C}-\text{CH}_3 \\ \text{I} \\ \text{OH} & \text{CH}_3 \\ \end{array}$$


 M_{r}


313.7

Concentration

2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	302 nm 247 nm		296 nm 243 nm	296 nm 242 nm
E ^{1%} 1cm	104 343		77 279	81 286
ε	3270 10760		2420 8750	2560 8960

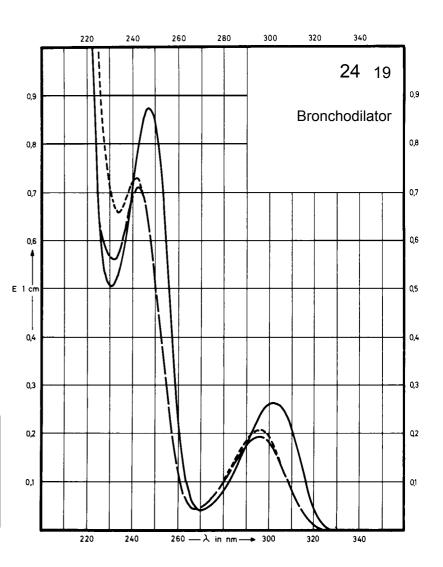

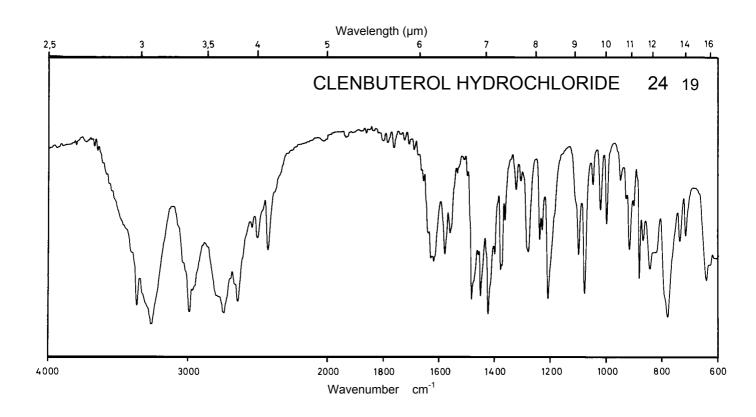



Name CAMYLOFINE DIHYDROCHLORIDE

 M_{r}


393.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm 262 nm 258 nm	268 nm 262 nm 257 nm	268 nm 262 nm 257 nm	
E ^{1%} 1cm	3.88 5.60 6.15	3.73 5.68 6.35	5.93 8.10 7.45	
ε	152 220	147 223	233 319	

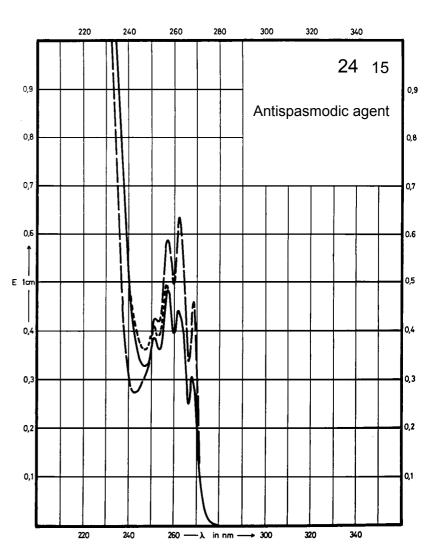

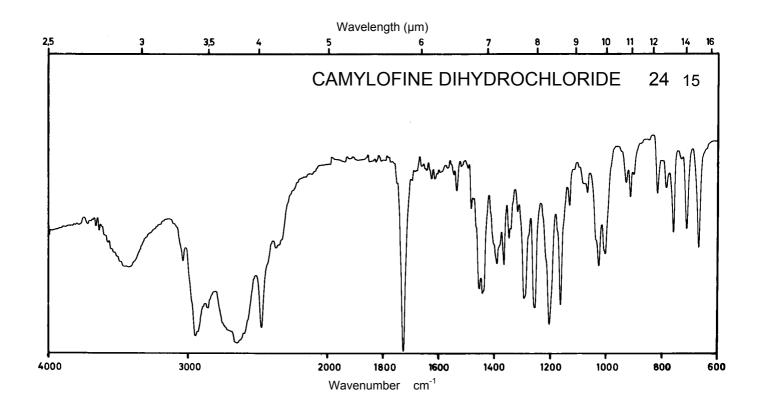

Name CHLORMEZANONE

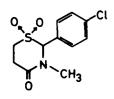



 M_{r}

237.7


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 259 nm 227 nm		Decom- position observed	Decom- position observed
E ^{1%} 1cm	10.4 9.3 637			
ε	247 220 15150			

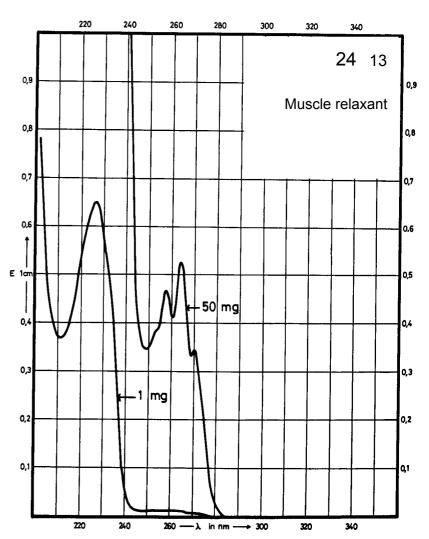

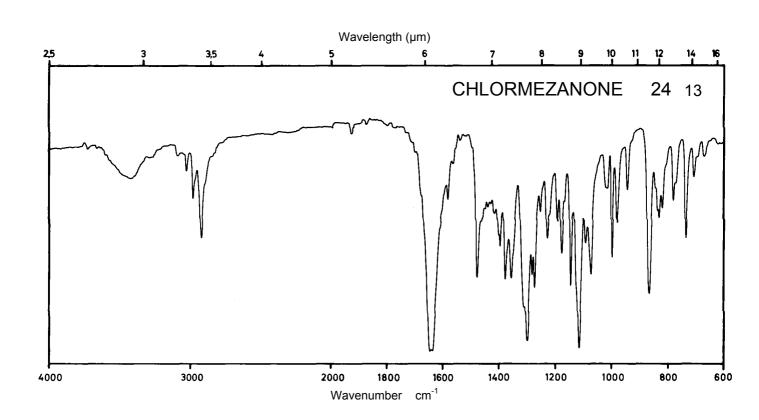



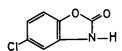



169.5

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm		278 nm	288 nm 243 nm
E ^{1%} 1cm	326		312	398 556
3	5530		5290	6750 9420

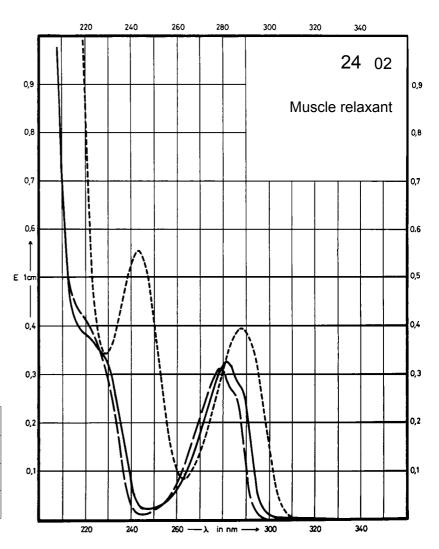

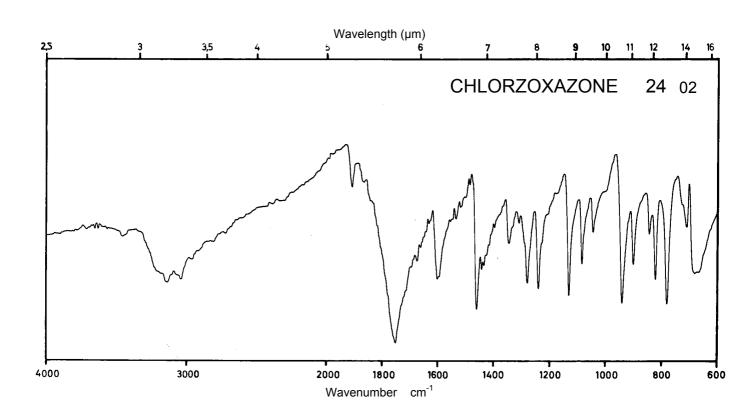


Name CAROVERINE HYDROCHLORIDE

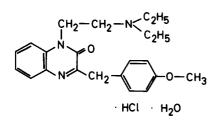



 M_{r}


419.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	332 nm 282 nm		334 nm 282 nm	340 nm 283 nm
E ^{1%} 1cm	184 197		192 183	190 178
3	7730 8270		8060 7680	7980 7470

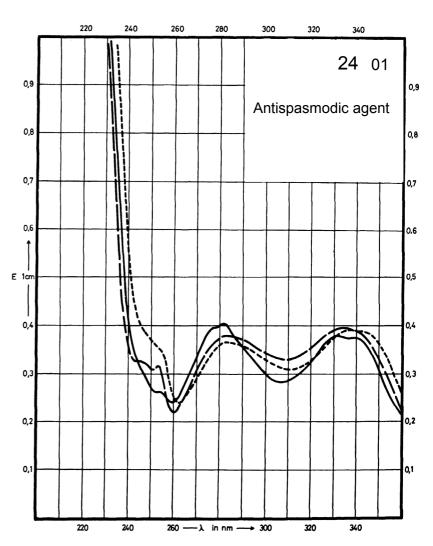

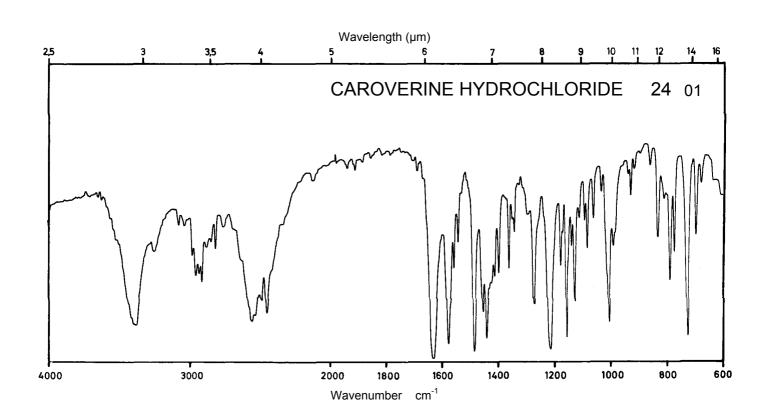

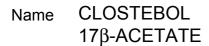

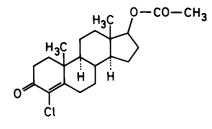



 M_r

Concentration 1.07 mg / 100 ml


364.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	254 nm			
E ^{1%} 1cm	374			
ε	13660			

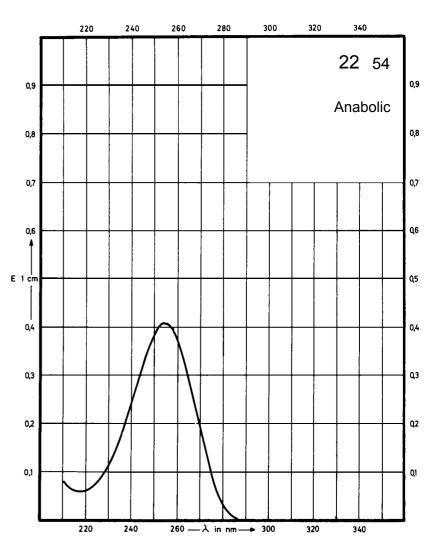

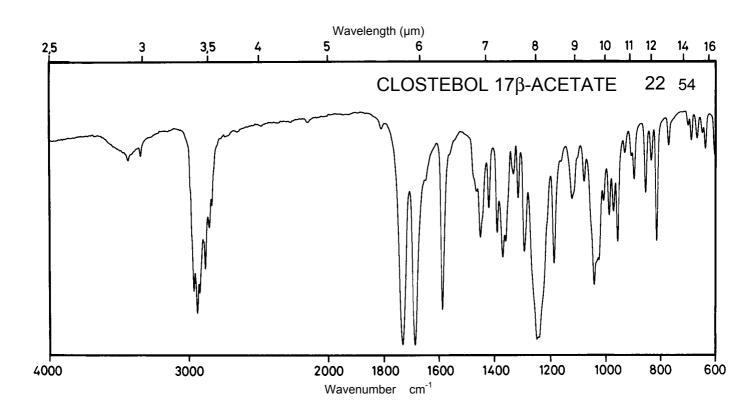





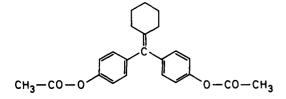


364.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	247 nm			
E ^{1%} 1cm	462			
ε	16830			

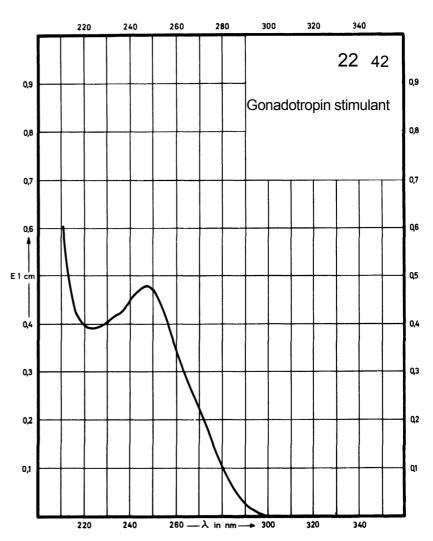

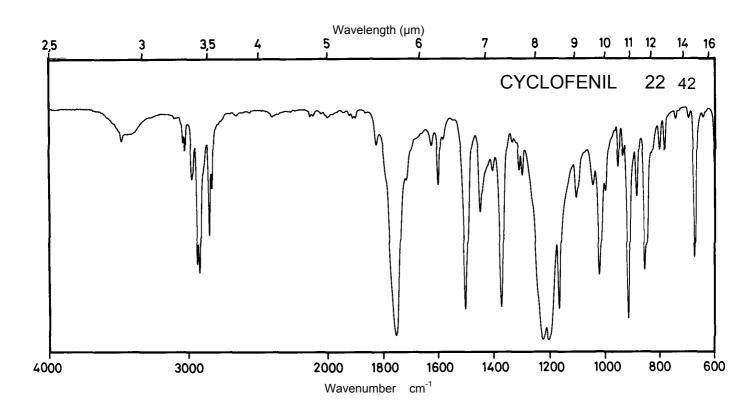





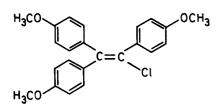


380.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	306 nm 246 nm			
E ^{1%} 1cm	412 651			
3	15700 24800			

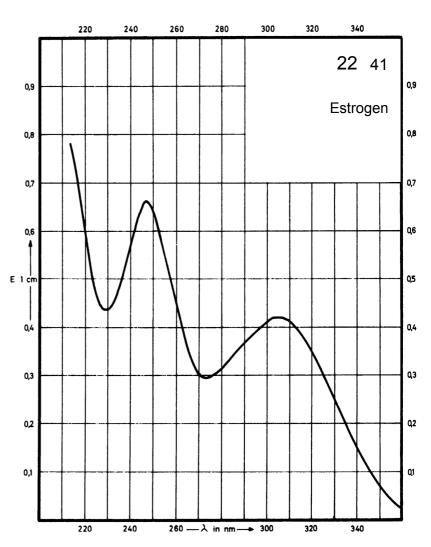

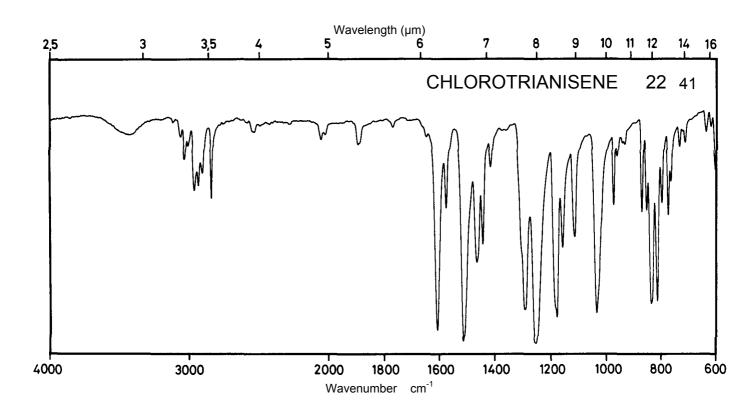




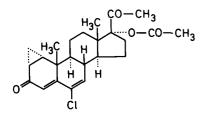


 M_r


416.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm	290 nm		
E ^{1%} 1cm	413	384		
3	17200	16000		

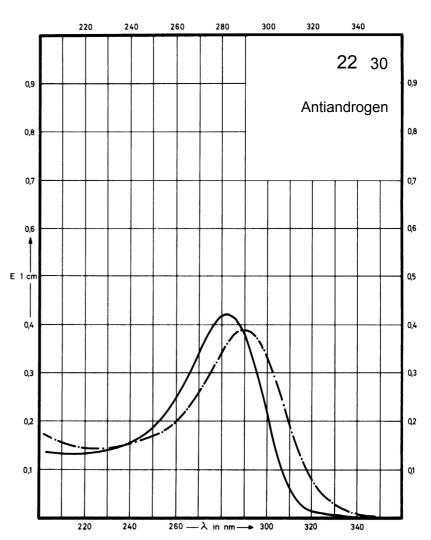

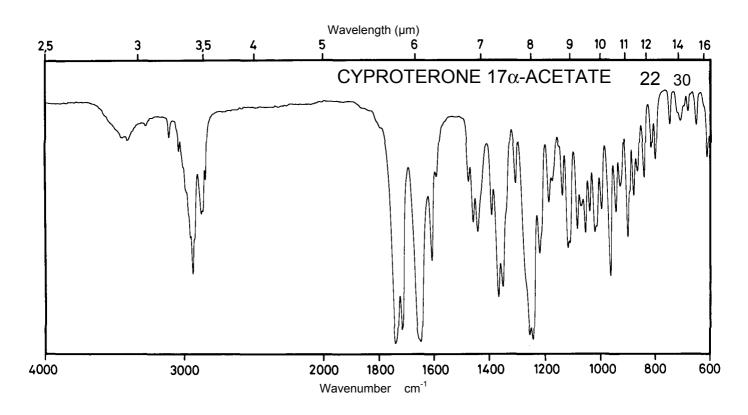





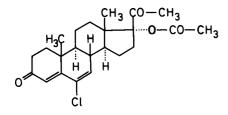


404.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm	292 nm		
E ^{1%} 1cm	548	527		
ε	22200	21300		

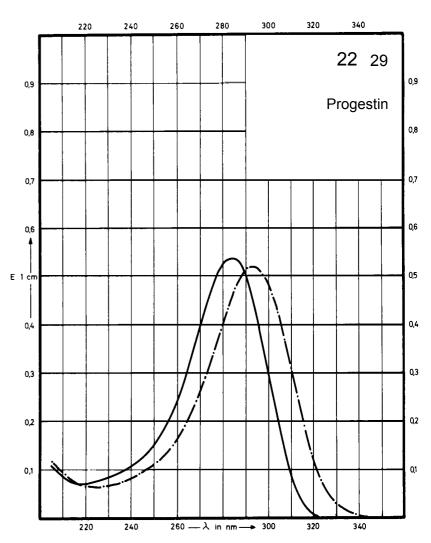

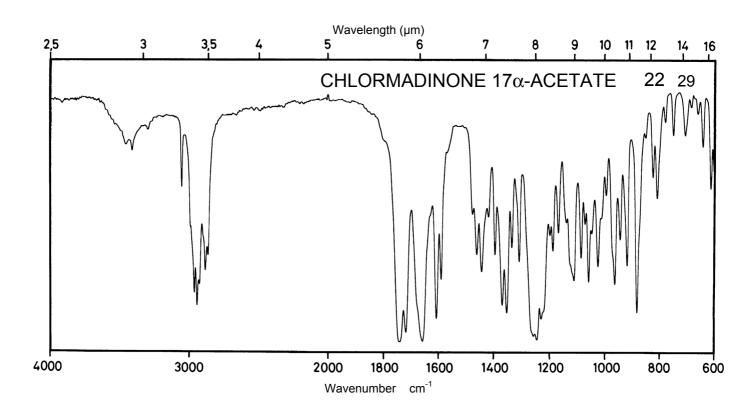



Name CYPROTERONE

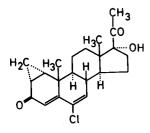


 M_r 374.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm	291 nm		
E ^{1%} 1cm	450	423		
3	16870	15840		

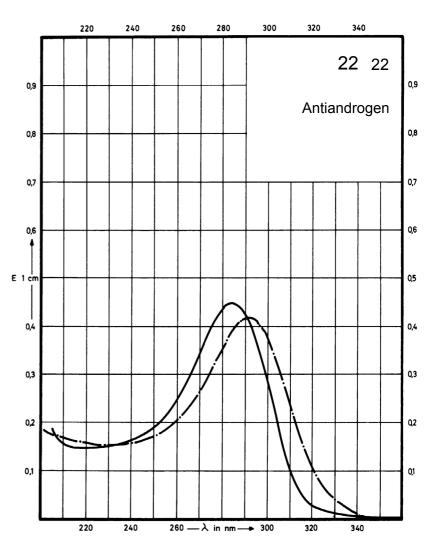

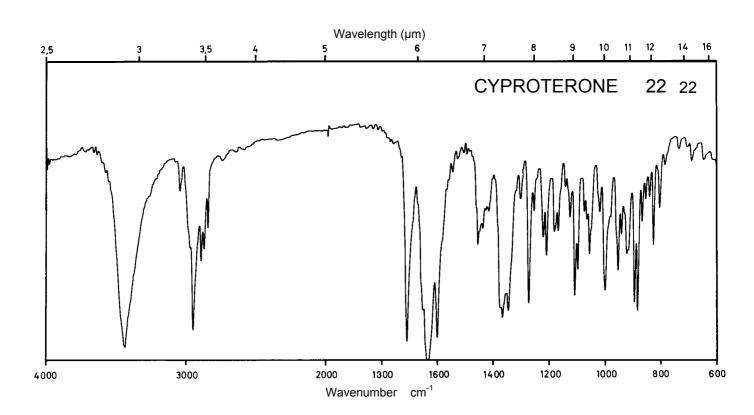



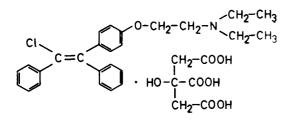



598.1

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	295 nm 234 nm	290 nm 234 nm	290 nm 234 nm	
E ^{1%} 1cm	193 325	175 312	175 312	
3	11540 19440	10470 18660	10470 18660	

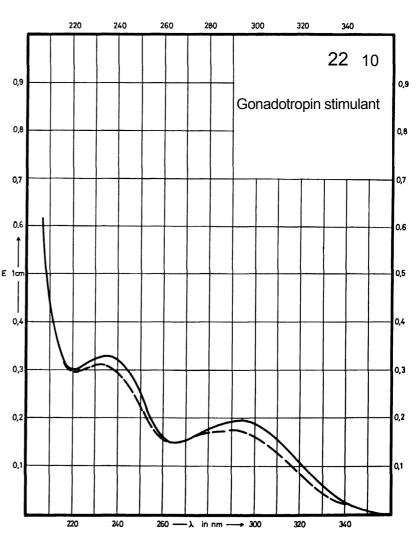

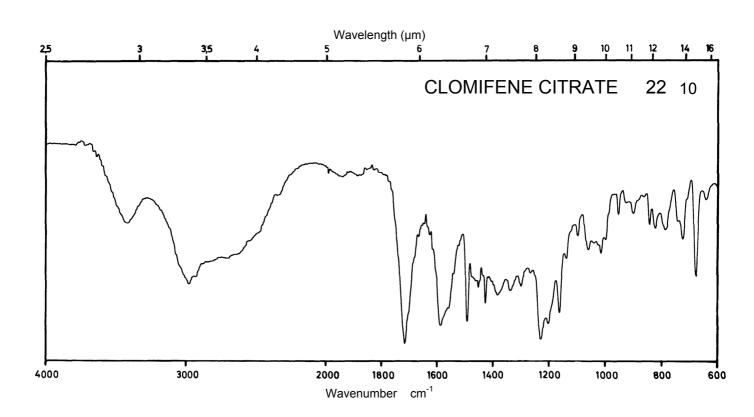





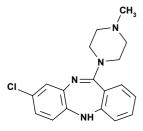


326.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	295 nm 259 nm 228 nm	293 nm 214 nm	295 nm 240 nm	293 nm
E ^{1%} 1cm	337 550 848	338 993	309 802	352
3	11000 18000 27700	11100 32400	10100 26200	11500

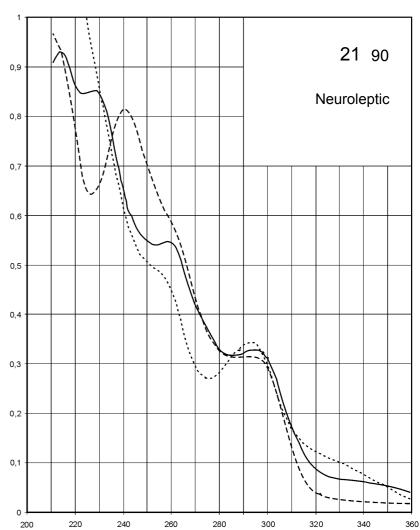

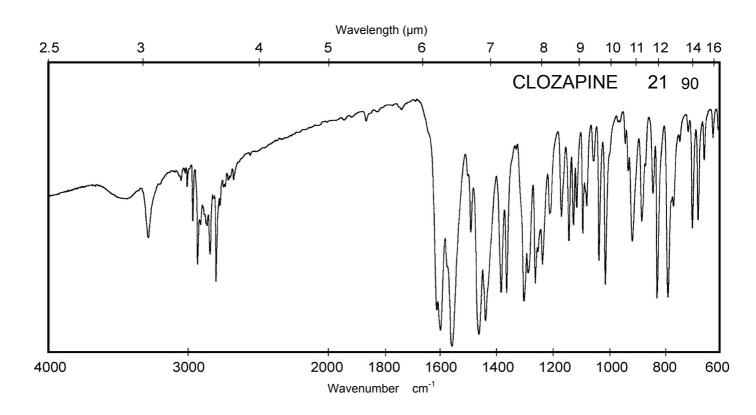




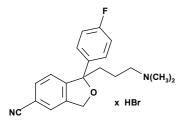


405.3

Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm 239 nm	239 nm	239 nm	239 nm
E ^{1%} 1cm	45 389	394	402	398
3	1820 15750	16000	16300	16100

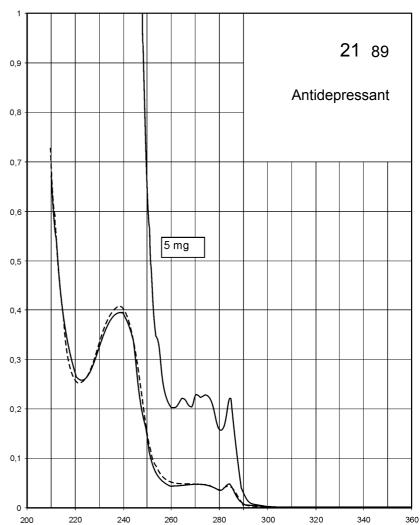


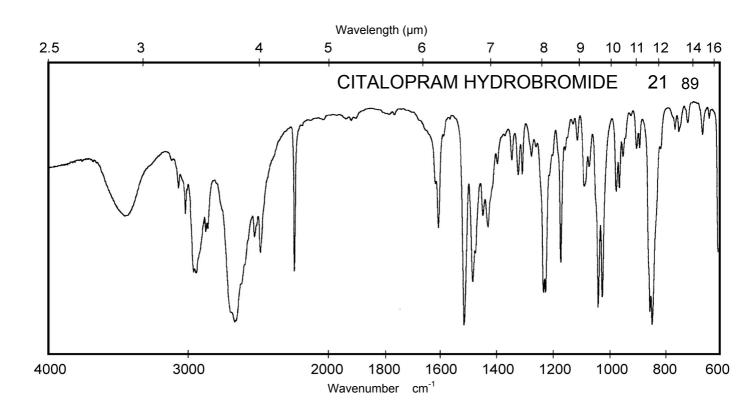


Name CLOPENTHIXOL DIHYDROCHLORIDE

 M_r 473.9

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	325 nm	322 nm	323 nm	322 nm
	270 nm	268 nm	268 nm	268 nm
absorption	230 nm	230 nm	230 nm	200 11111
1%	64	59	61	55
E ^{1%} 1cm	297	282	287	258
TCM	707	692	700	230
	3040	2800	2900	2600
3	14000	13400	13600	12200
	33500	32800	33200	12200

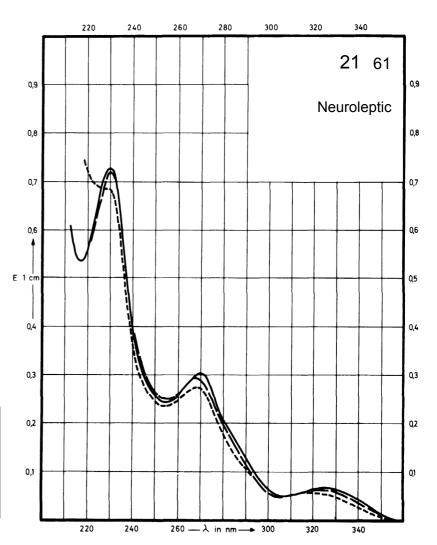

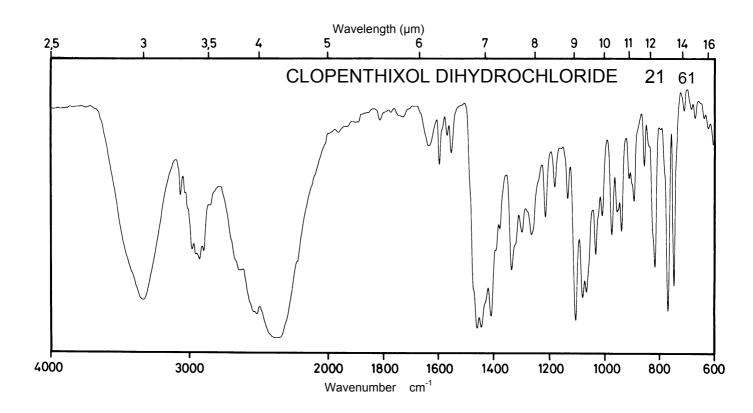



Name CLOMETHIAZOLE EDISILATE

$$\begin{bmatrix} C_1 - CH_2 - CH_2 & S \\ CH_3 & N \end{bmatrix}_2 \cdot C_2H_6O_6S_2$$


 M_r 513.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	250 nm	249 nm	257 nm	249 nm
E ^{1%} 1cm	175	174	184	170
3	8980	8910	9470	8730

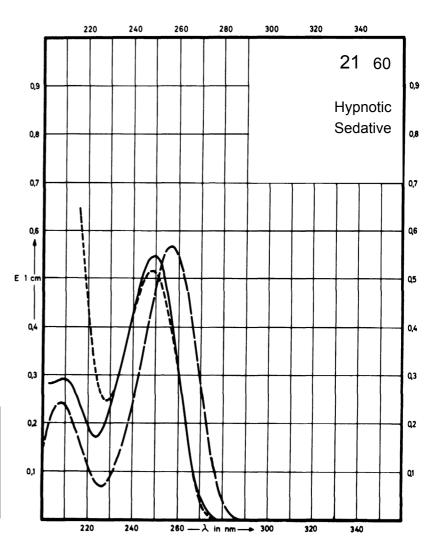

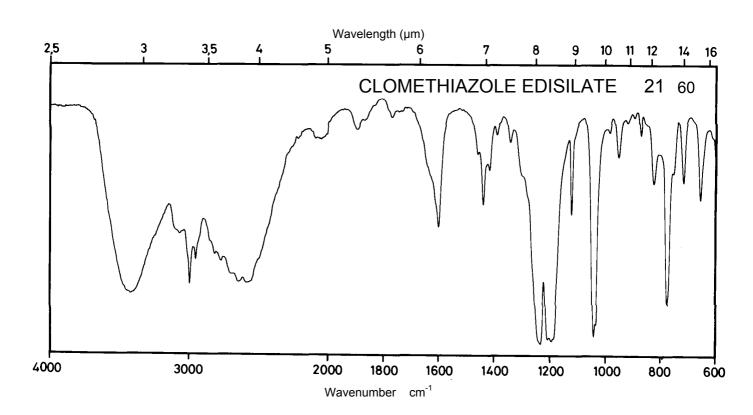




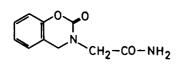


 M_r 206.2


Concentration 1 mg / 100 ml 2 mg / 100 ml 10 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Decom-position observed Maximum of 275 nm 268 nm 274 nm 267 nm 274 nm 267 nm absorption E^{1%} 1cm 35 39 41 44 35 39 845 910 720 805 ε

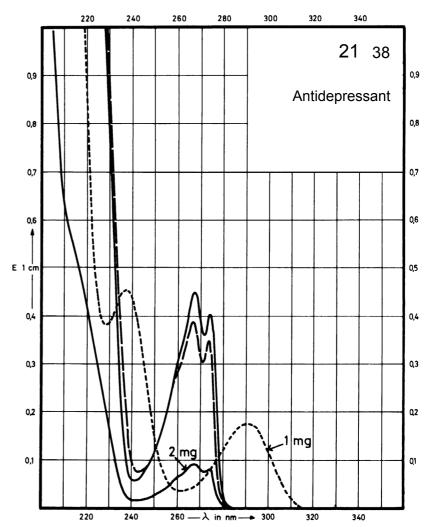

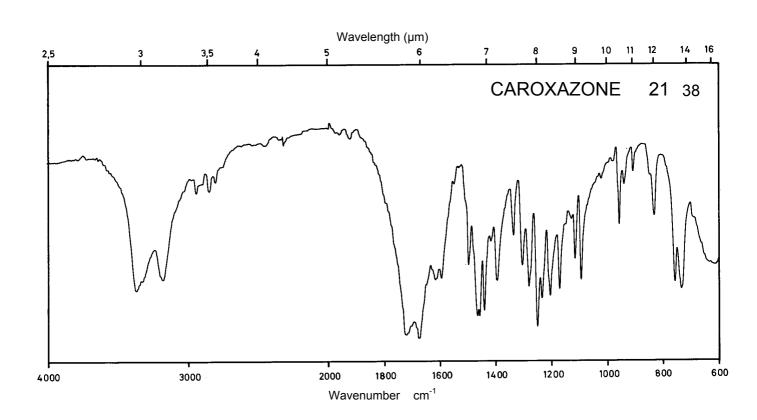


Name NORTRIPTYLINE HYDROCHLORIDE

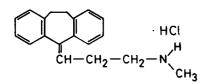


 M_{r}


299.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm		239 nm	
E ^{1%} 1cm	472		456	
ε	14160		13680	

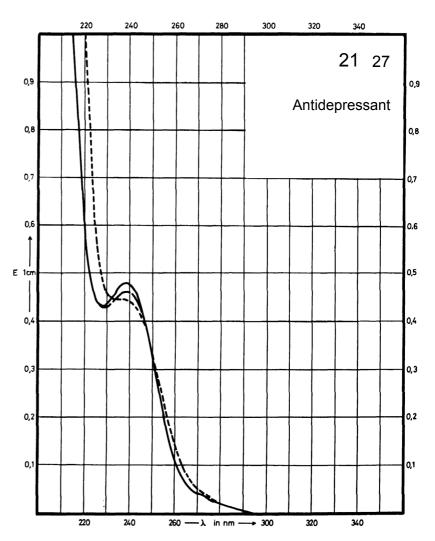

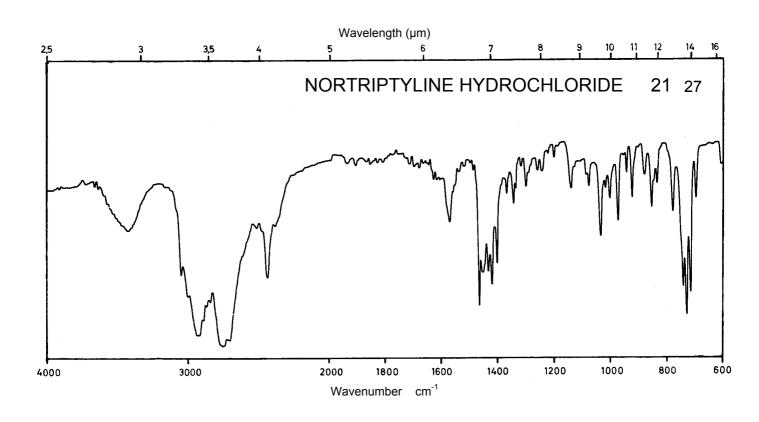


Name CLOMIPRAMINE HYDROCHLORIDE

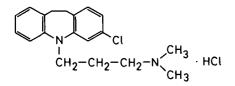


 M_{r}

351.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	251 nm		251 nm	273 nm 254 nm
E ^{1%} 1cm	238		228	220 237
ε	8360		8020	7730 8330

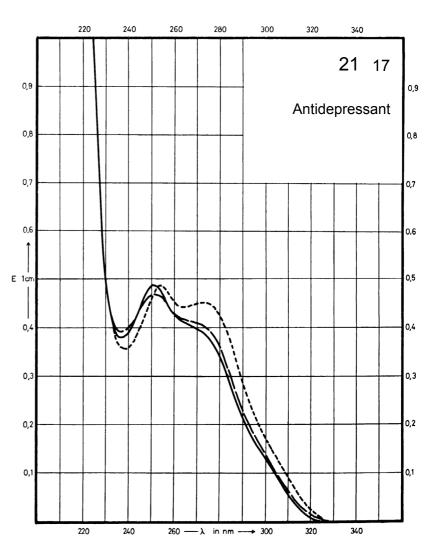

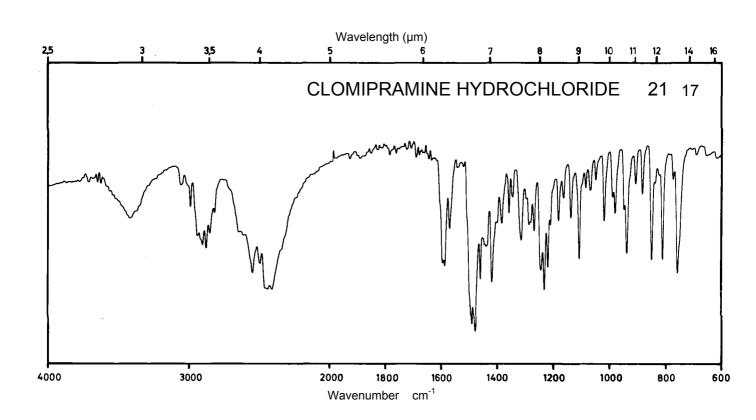



Name CLOFENCICLAN HYDROCHLORIDE

 M_{r}


346.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm 266 nm 274 nm		259 nm 265 nm	
E ^{1%} 1cm	5.84 6.65 4.00		6.19 6.94	
3	202 230 139		214 240	

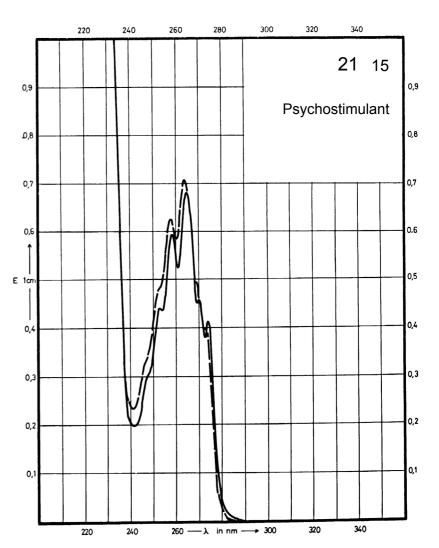

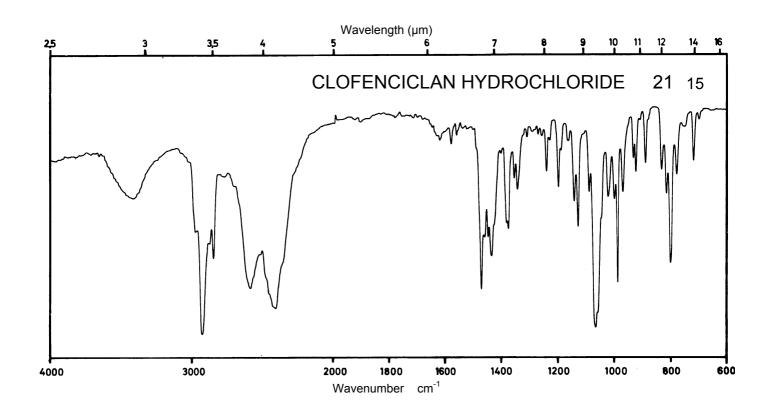



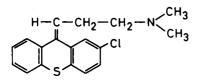



 M_{r}

315.9


Concentration 0.5 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	326 nm 270 nm		324 nm 268 nm 228 nm	
E ^{1%} 1cm	100 449		98 443 1096	
ε	3160 14180		3100 13990 34620	

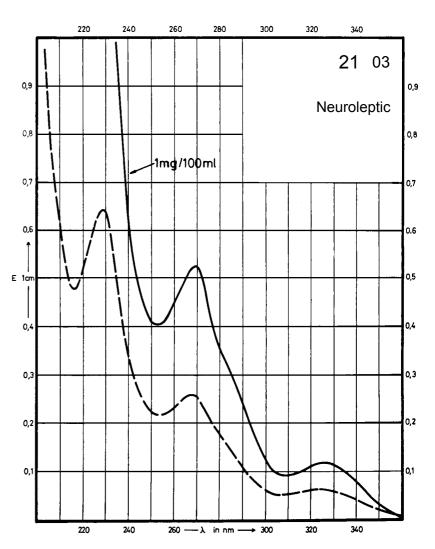

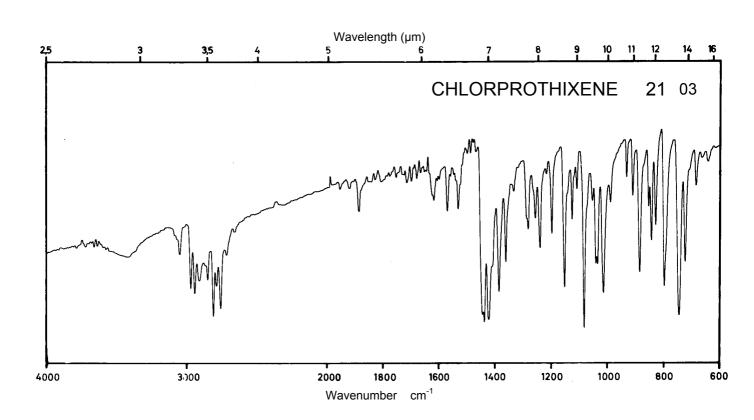

Name CHLORPROMAZINE HYDROCHLORIDE

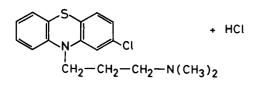


 M_{r}


355.4

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm 255 nm		306 nm 254 nm	
E ^{1%} 1cm	118 1002		115 950	
ε	4190 35610		4090 33760	

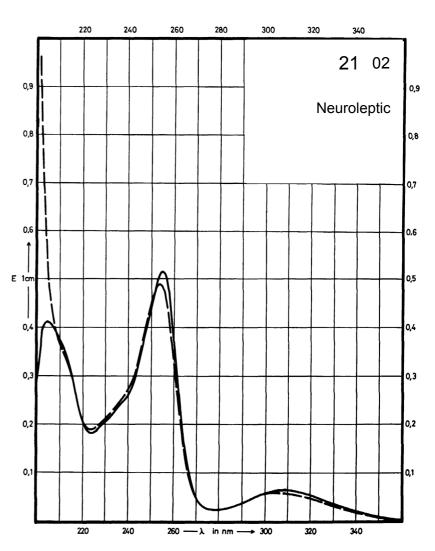

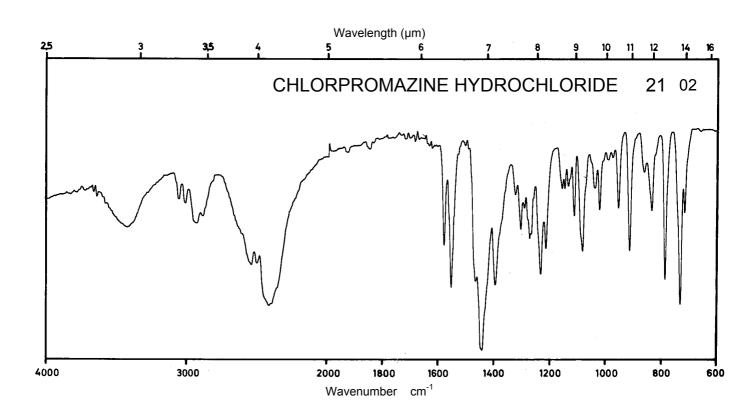


Name CELIPROLOL HYDROCHLORIDE

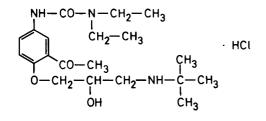


 M_r


416.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	330 nm 233 nm		326 nm 231 nm	326 nm 230 nm
E ^{1%} 1cm	60 786		59 668	62 660
ε	2500 32700		2470 27800	2580 27500

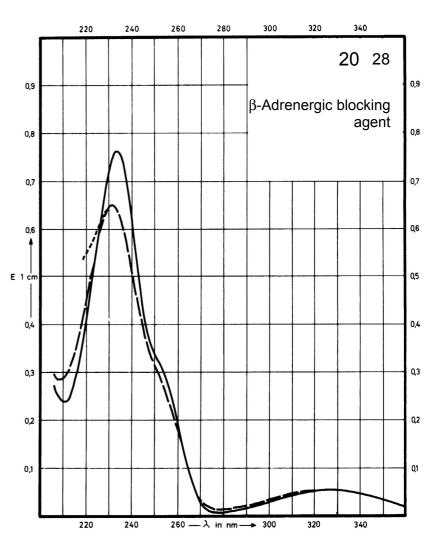

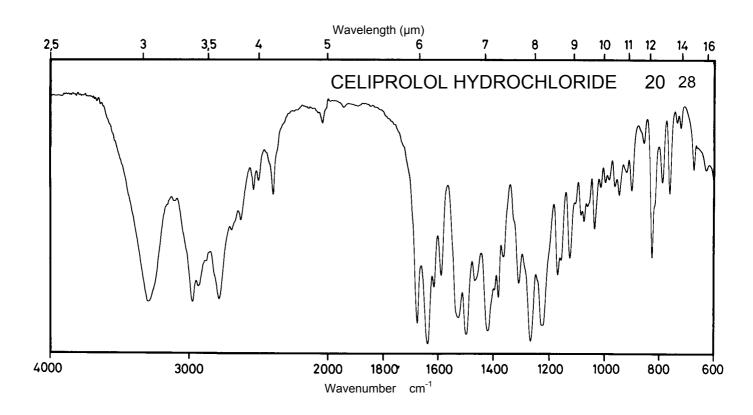



Name CARTEOLOL HYDROCHLORIDE

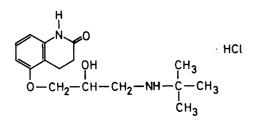


 M_r 328.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 279 nm 251 nm		251 nm	251 nm
E ^{1%} 1cm	86 91 308		300	295
ε	2800 3000 10100		9900	9700

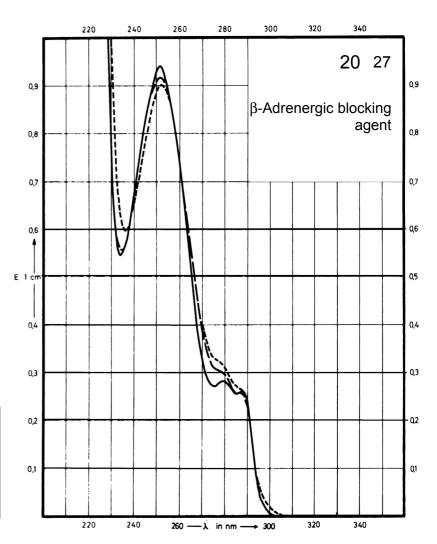

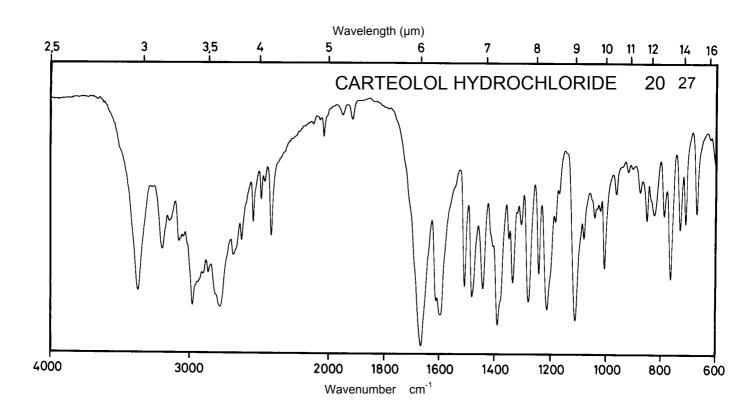



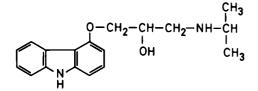


 M_r


298.4

Concentration

0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	332 nm 285 nm 242 nm	331 nm 284 nm 240 nm	331 nm 284 nm 240 nm	331 nm 284 nm 240 nm
E ^{1%} 1cm	209 495 1760	170 388 1530	177 410 1620	180 410 1600
ε	6240 14770 52500	5070 11580 45660	5280 12230 48340	5370 12230 47740

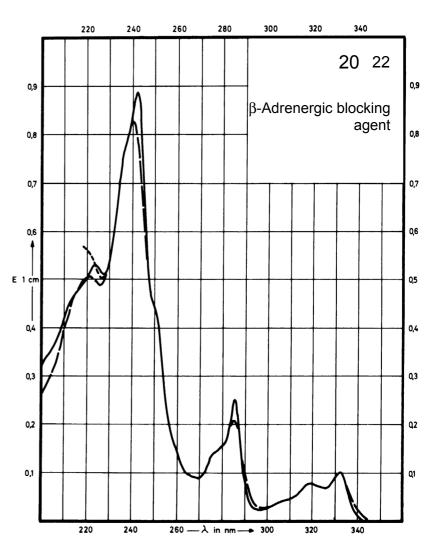

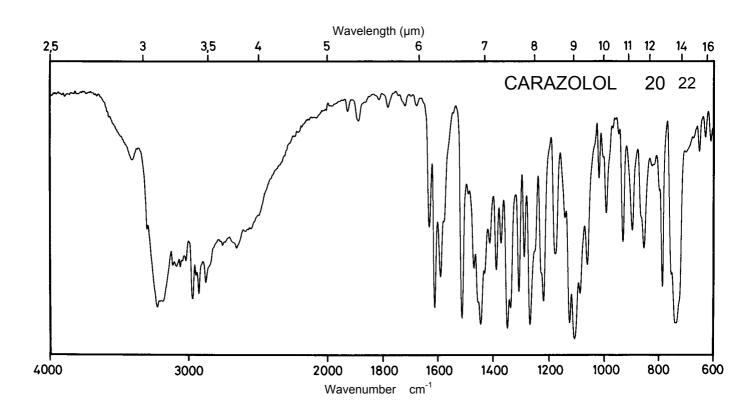


Name CLIDINIUMBROMID


Mol.-Gew.

432,4


Gemessene Konz. 53 mg / 100 ml

Lösungsmittel Symbol	Methanol	Wasser	0,1 M-HCI	0,1 M-NaOH
Absorptions- maximum	265 nm 259 nm 253 nm	259 nm 253 nm	259 nm 253 nm	
E ^{1%} _{1cm}	8,28 10,6 9,34	11,3 10,1	11,3 10,1	
ε	358 457 412	490 437	490 437	

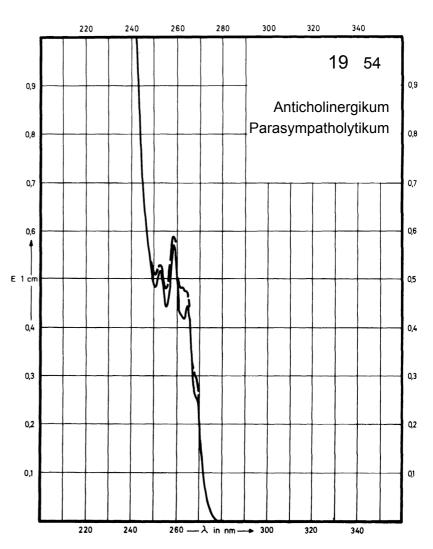

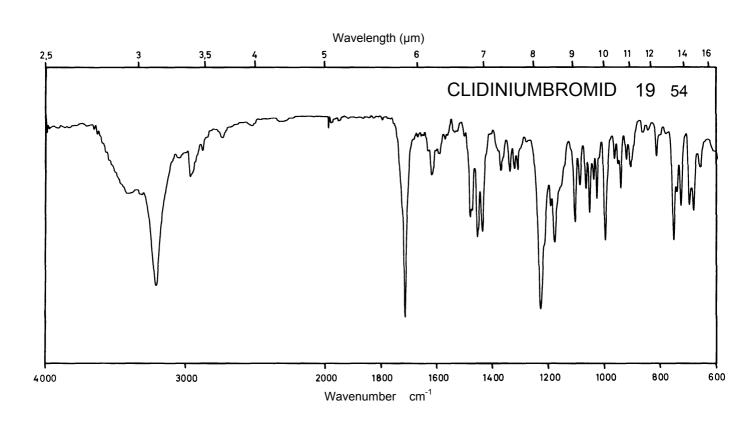




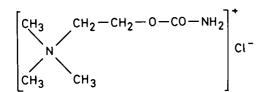


 M_{r}


182.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

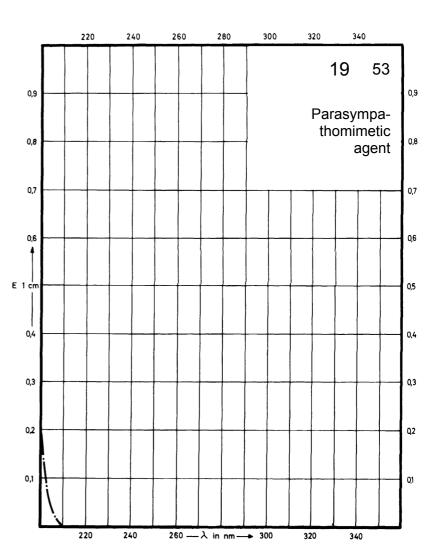

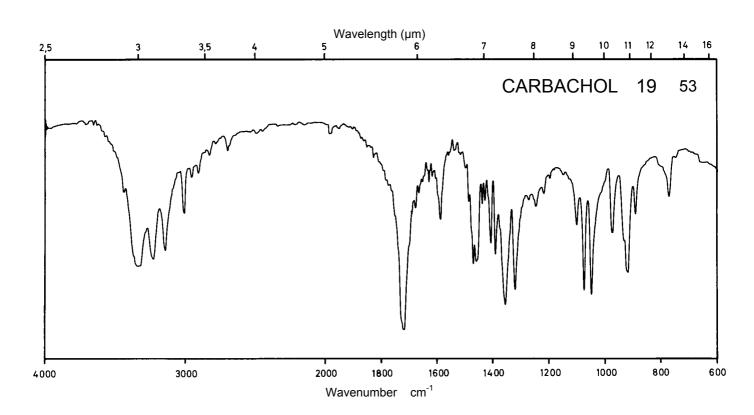





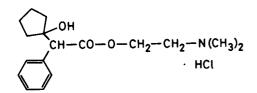


 M_r 327.9

Concentration 5 mg / 100 ml 100 mg / 100 ml

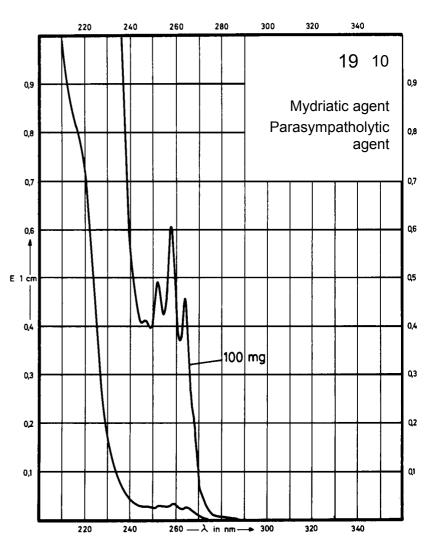

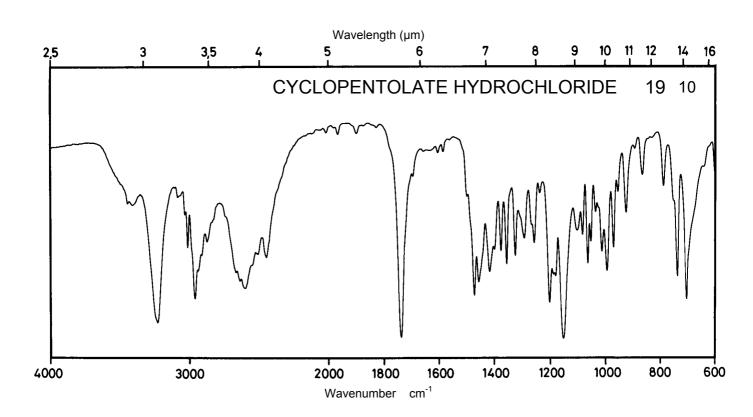




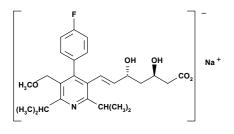


 M_{r}


481.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		267 nm	285 nm	
E ^{1%} 1cm		118	196	
ε		5670	9440	

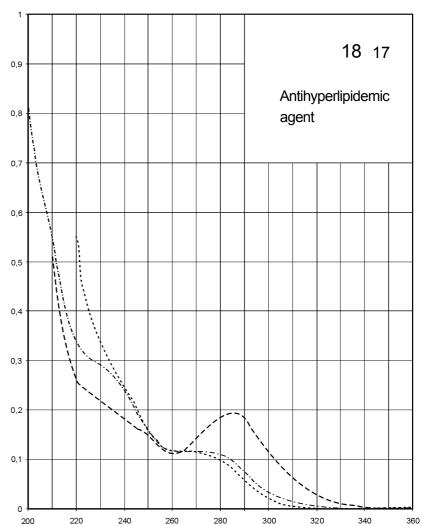

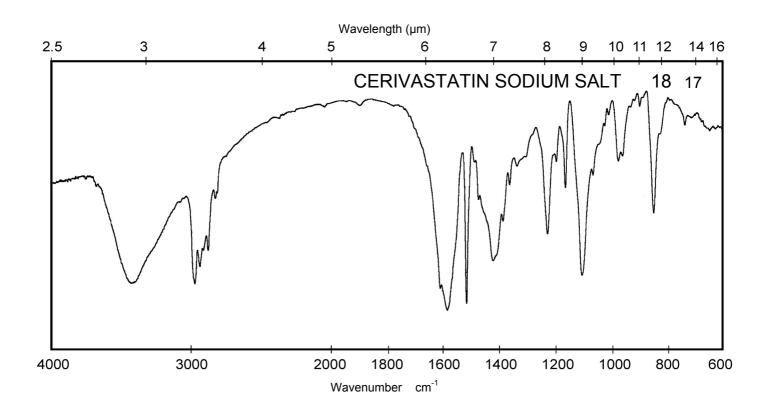




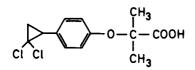


 M_r 289.2


Concentration 2 mg / 100 ml 15 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm	275 nm	273 nm	275 nm
	232 nm	233 nm	230 nm	233 nm
E ^{1%} 1cm	30	33	26	36
	441	446	425	463
ε	870	960	740	1000
	12750	12900	12300	13400

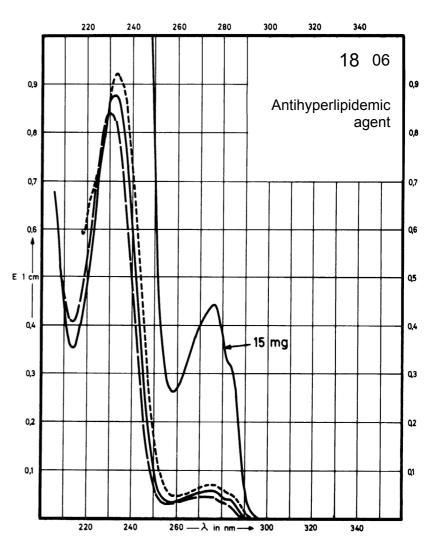

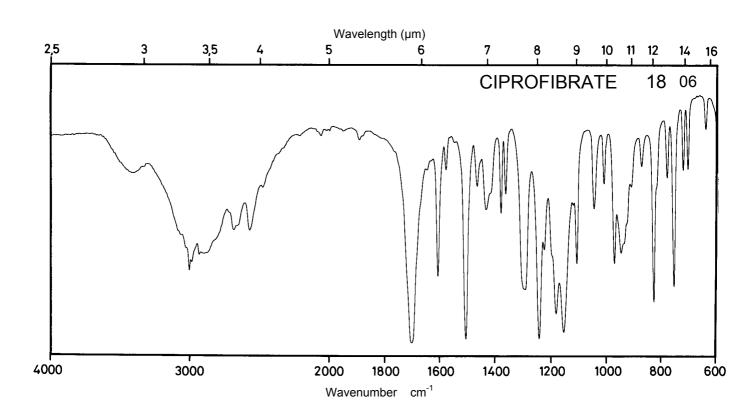




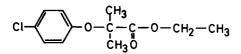


 M_r 242.7

Concentration 0.9 mg / 100 ml 9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	279 nm 226 nm			
E ^{1%} 1cm	42 463			
3	1020 11240			

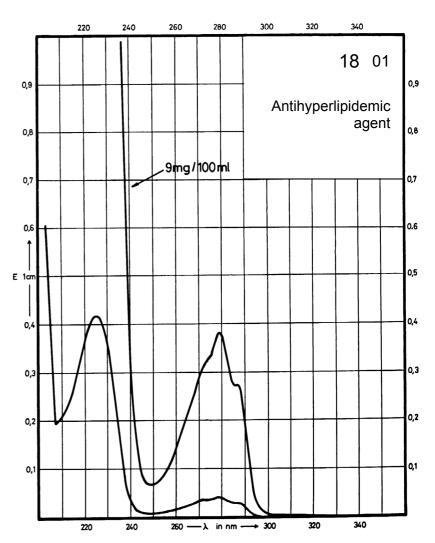

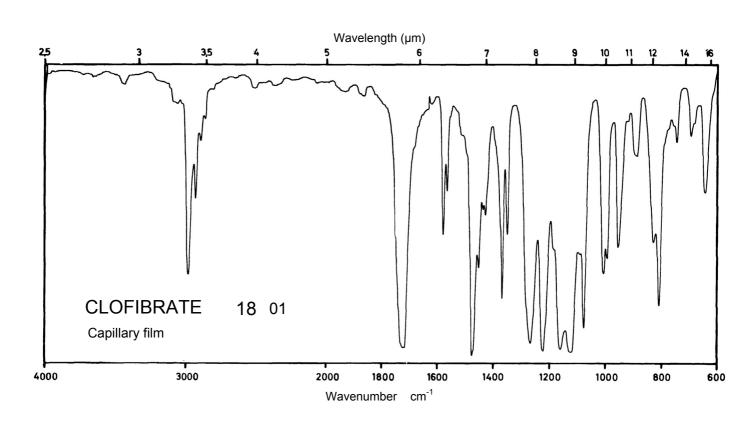


Name COCAINE HYDROCHLORIDE


 $\,M_r\,$

339.8


Concentration 2 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm 230 nm	275 nm 232 nm	275 nm 232 nm	269 nm
E ^{1%} 1cm	28.4 405	30 390	30 390	21
ε	960 13740	1020 13240	1020 13240	710

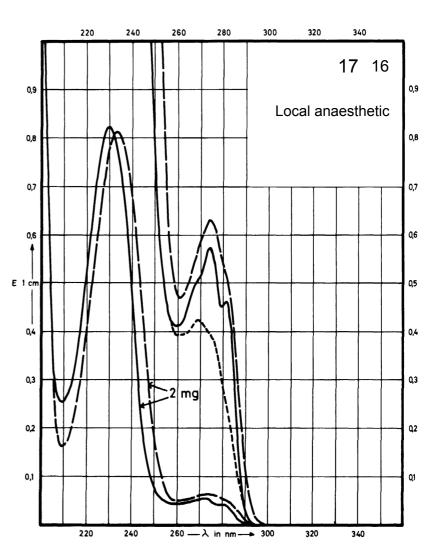

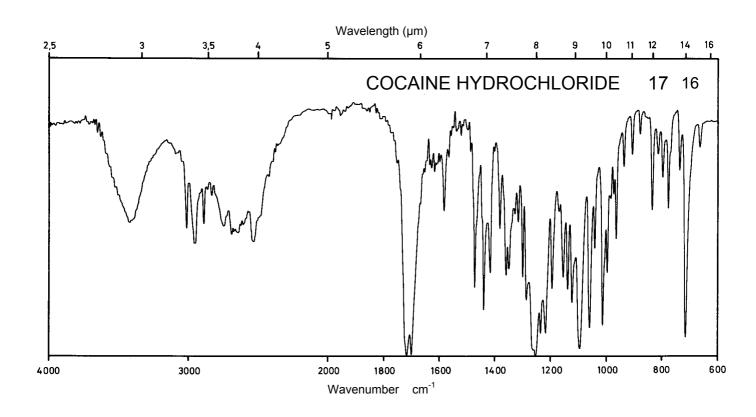

Name CINCHOCAINE HYDROCHLORIDE

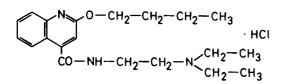



 M_r


379.9


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	328 nm		318 nm 247 nm	328 nm 239 nm
E ^{1%} 1cm	111		215 625	111 485
3	4220		8180 23750	4220 18420

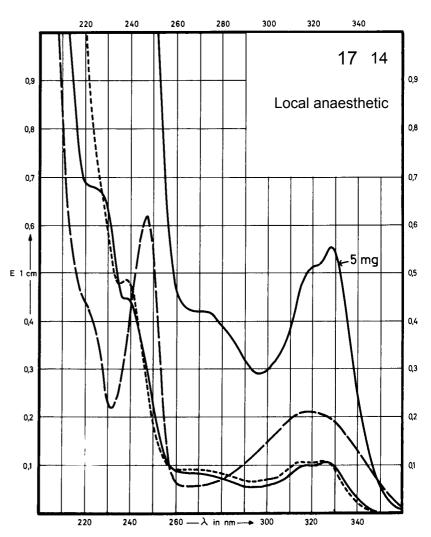

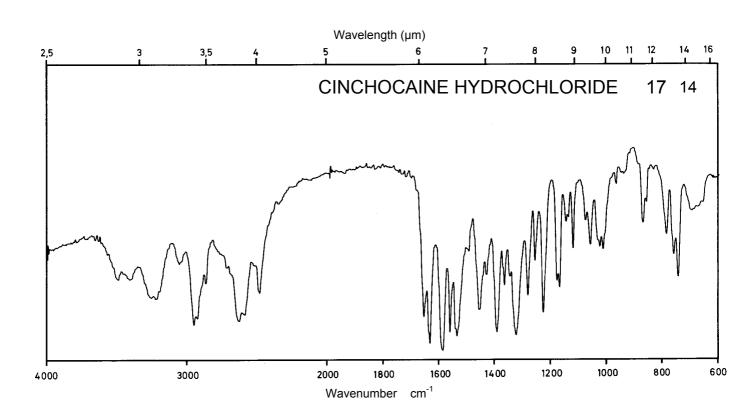

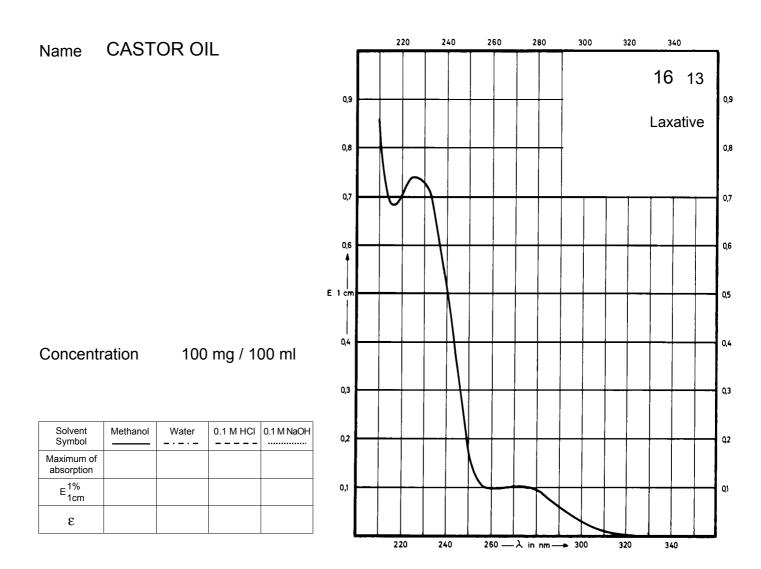

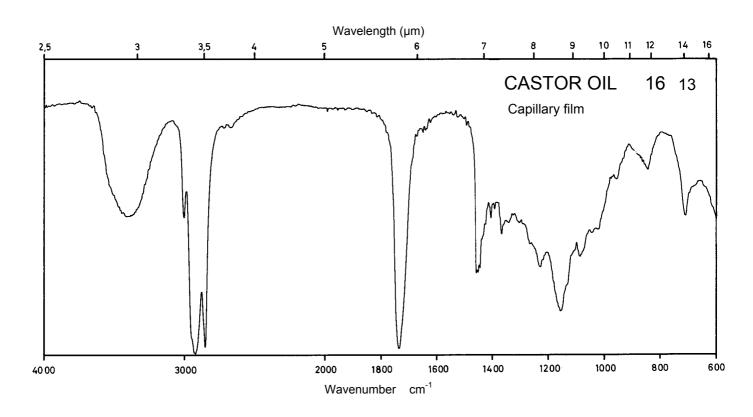



Name CIMETIDINE

$$H_3C$$
 H_3C
 H_3C


 M_r 252.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	219 nm		219 nm	
E ^{1%} 1cm	810		790	
ε	20440		19930	

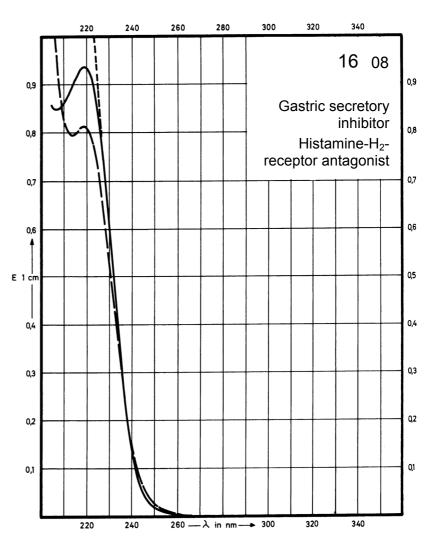

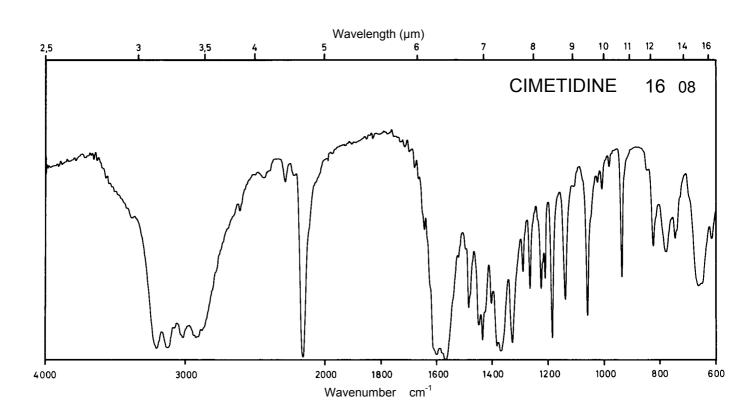




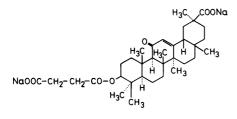


 M_r


614.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	252 nm			260 nm
E ^{1%} 1cm	172			172
ε	10600			10600

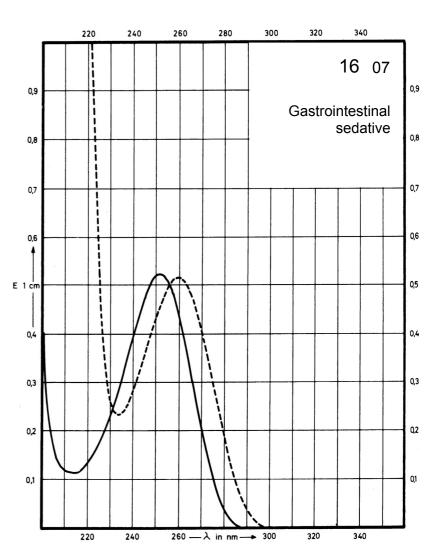

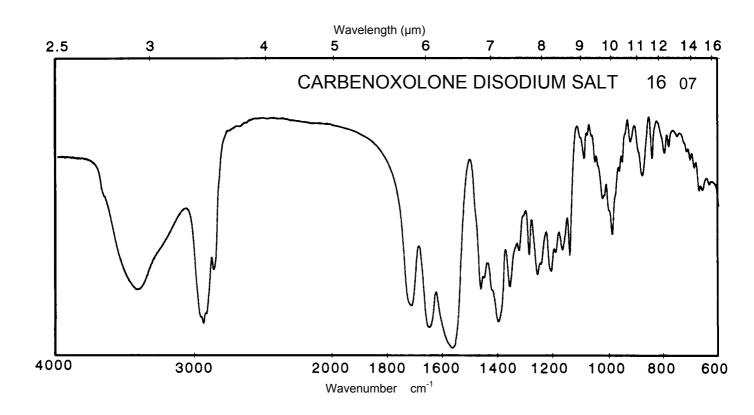





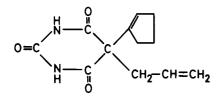

 M_r 234.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				244 nm
E ^{1%} 1cm				319
ε				2480

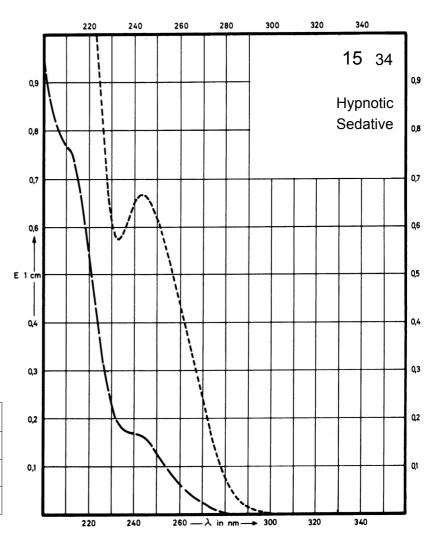

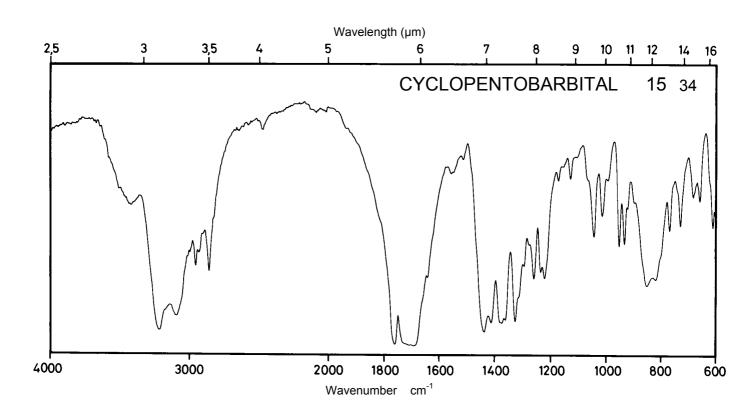


Name CARBROMAL


$$\begin{array}{c} \text{CH}_3 \\ \text{CH}_2 \\ \text{Br} - \text{C} - \text{CO} - \text{NH} - \text{CO} - \text{NH}_2 \\ \text{CH}_2 \\ \text{CH}_3 \end{array}$$

 M_r 237.1


Concentration 5 mg / 100 ml 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

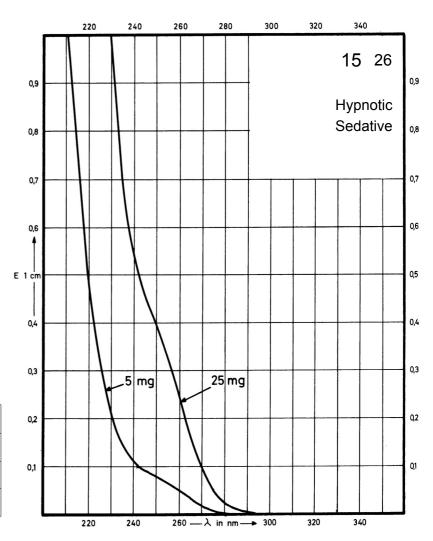

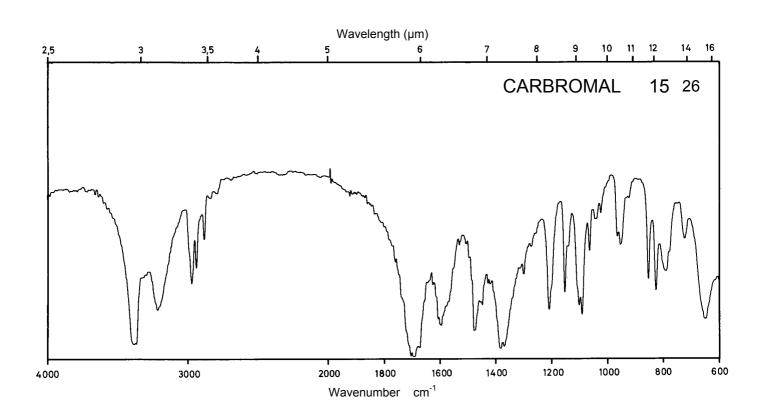


Name CLONAZEPAM

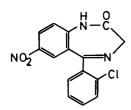


 M_{r}

315.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	309 nm 245 nm		274 nm	364 nm
E ^{1%} 1cm	361 460		645	448
3	11400 14540		20360	14130

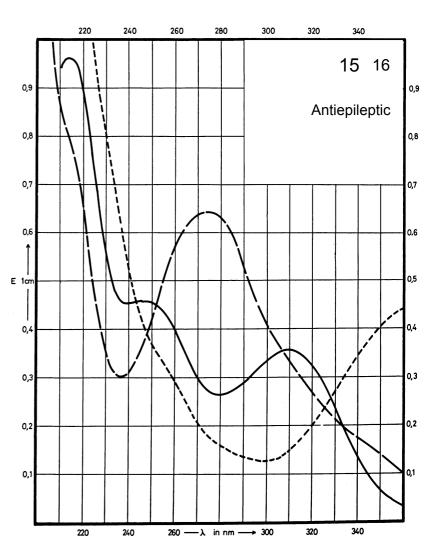

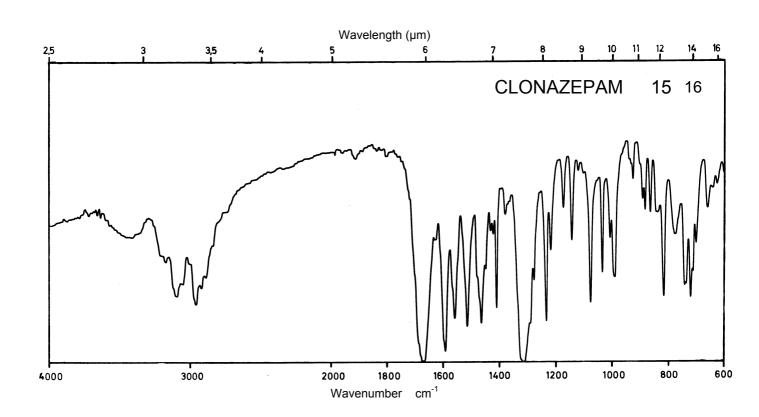

Name CARBAMAZEPINE

CO-NH₂


 M_{r}


236.3

Concentration

1.2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 237 nm		285 nm	285 nm
E ^{1%} 1cm	454 568		480	480
ε	10740 13410		11330	11330

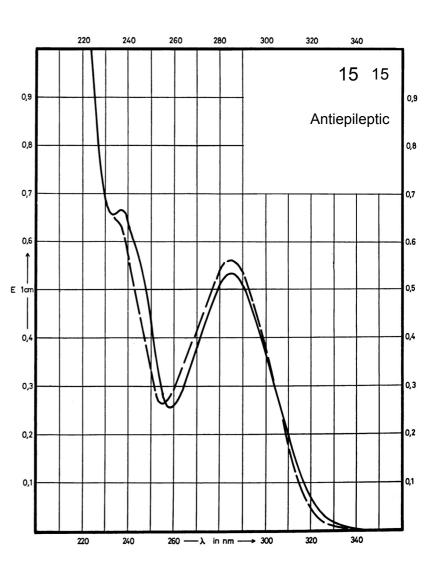

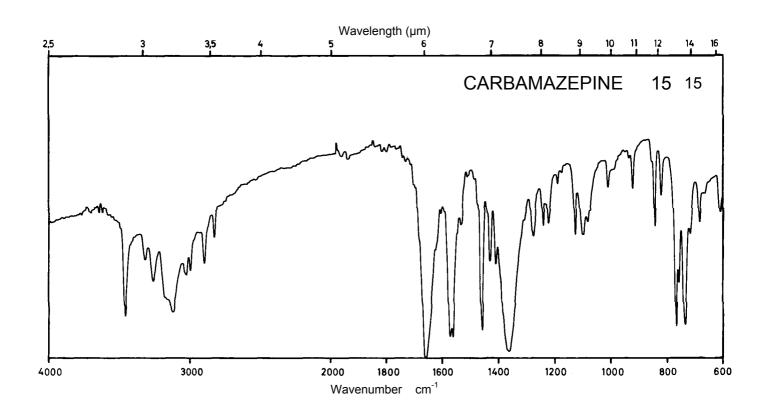


Name CYCLOBARBITAL CALCIUM SALT


$$\begin{bmatrix} & & & & & \\ & & & & \\ & & & & \\ & & & & \\ & & & & \\ & & & & \\ & & & & \\ & & \\ & & & \\ & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & &$$

 M_{r}


510.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				250 nm
E ^{1%} 1cm				268
ε				13680

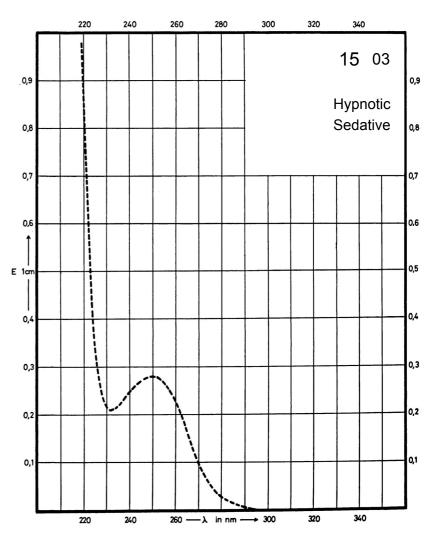

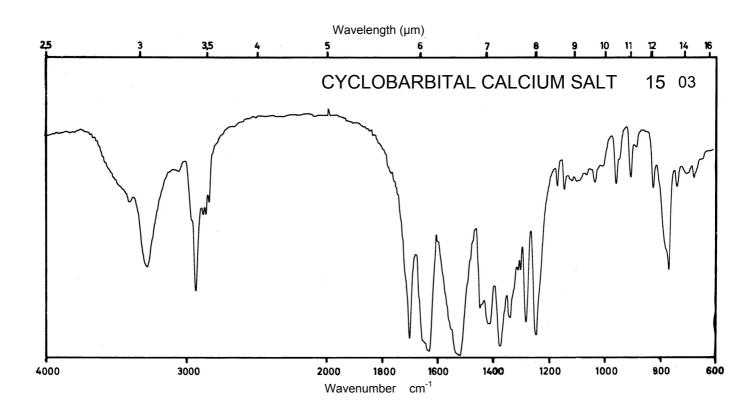



Name CARVEDILOL

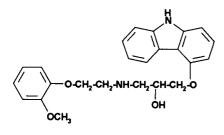


 M_r 406.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	332 nm	332 nm	332 nm	332 nm
	286 nm	285 nm	285 nm	285 nm
	242 nm	241 nm	241 nm	241 nm
E ^{1%} _{1cm}	163	116	121	123
	401	276	287	290
	1415	1030	1085	1065
3	6600	4700	4900	5000
	16300	11200	11700	11800
	57500	41800	44100	43300

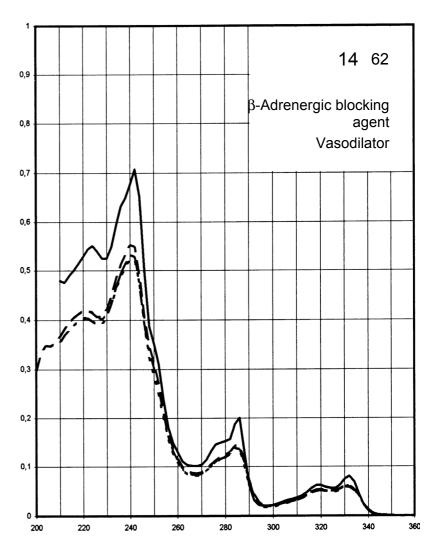

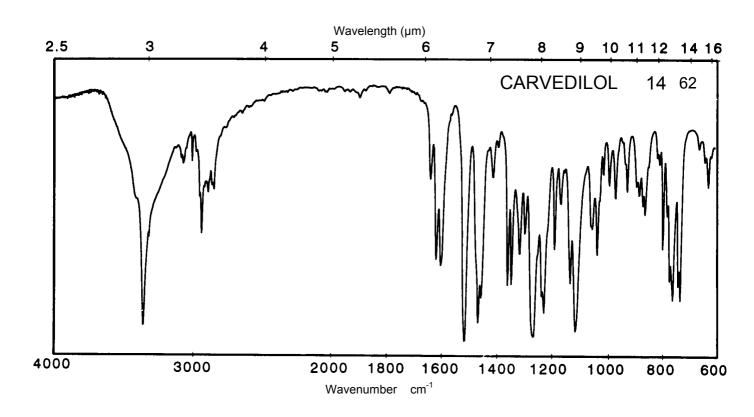



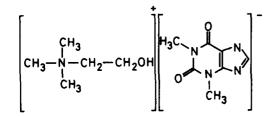


 M_{r}


283.3

Concentration 1.


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm		270 nm	275 nm
E ^{1%} 1cm	347		351	418
ε	9830		9930	11850

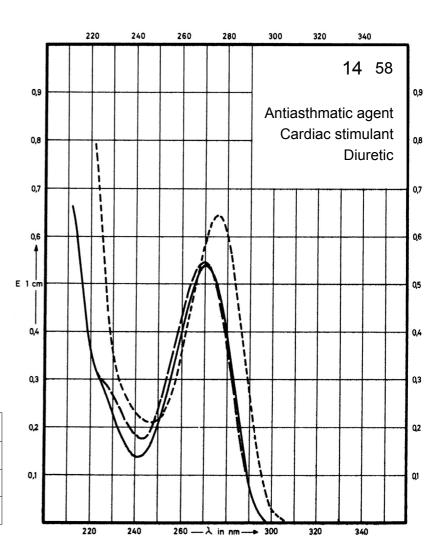

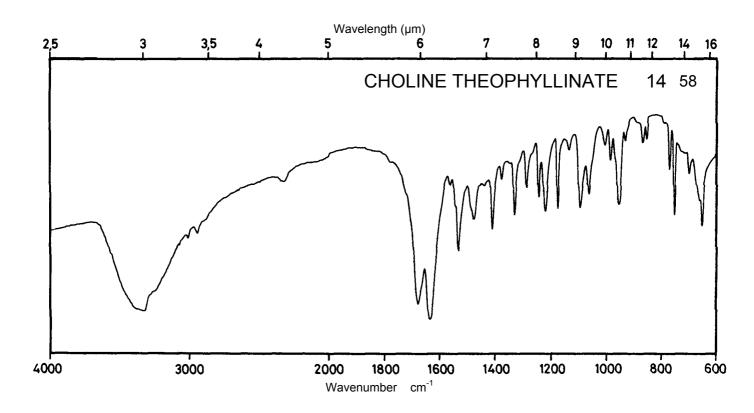




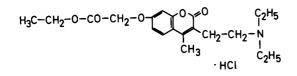


 M_{r}


397.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	321 nm		321 nm	
E ^{1%} 1cm	440		440	
3	17510		17510	

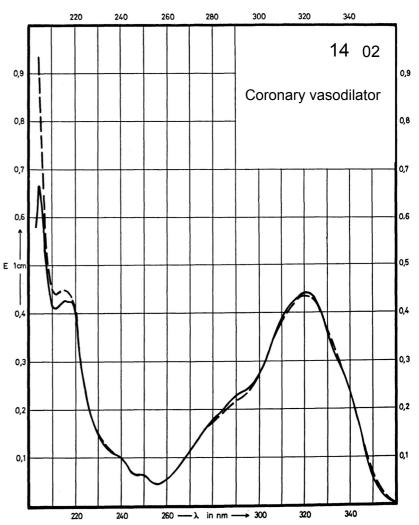

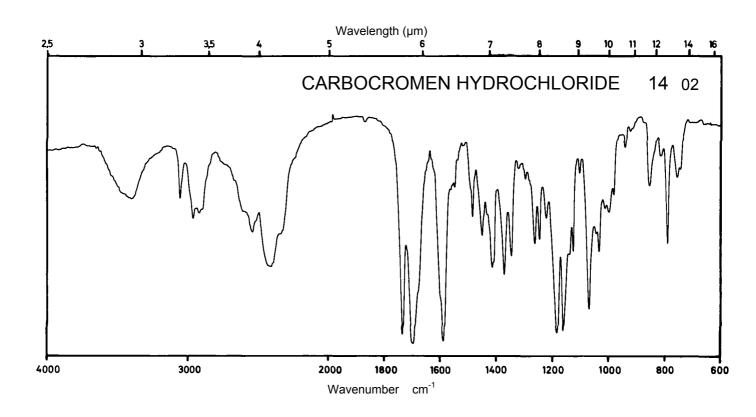



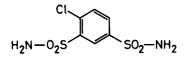

 M_r 270.7


Concentration 1 mg / 100 ml 10 mg / 100 ml


 Solvent Symbol
 Methanol —————
 Water —————
 0.1 M HCI —————
 0.1 M NaOH

 Maximum of absorption
 284 nm 275 nm 232 nm
 284 nm 276 nm 232 nm

 E 1% 1cm
 31 45 455
 41 45 455

 840 1100 950 12300
 1100 12300

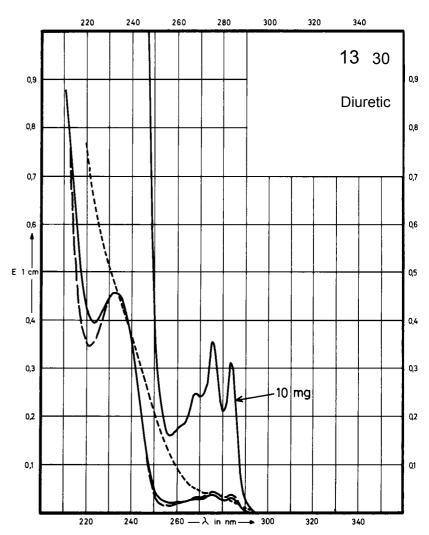

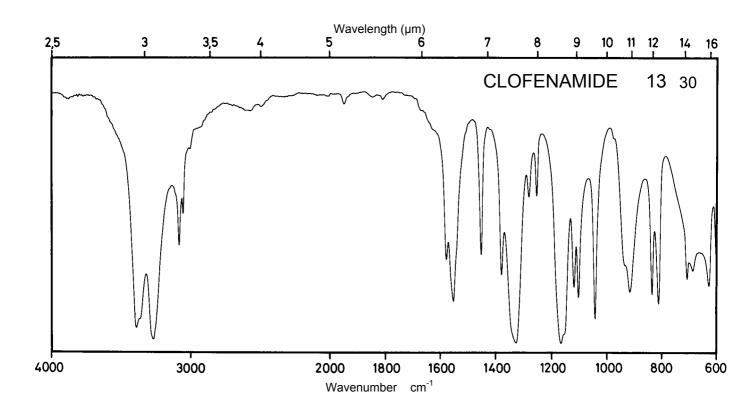


Name CYCLOPENTHIAZIDE


 M_{r}

379.9


Concentration 1 mg / 100 ml 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 272 nm		315 nm 272 nm	320 nm 273 nm
E ^{1%} 1cm	80 585		82 560	70 440
ε	3040 22220		3120 21270	2660 16720

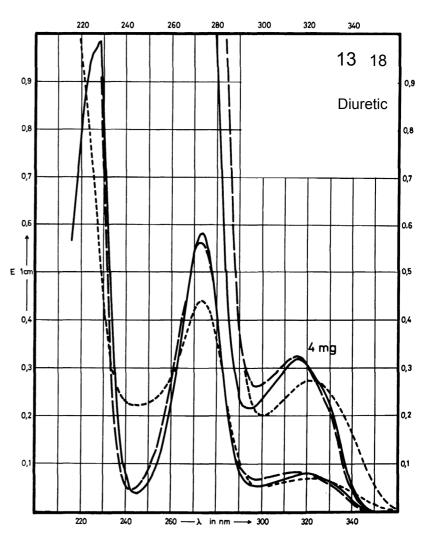

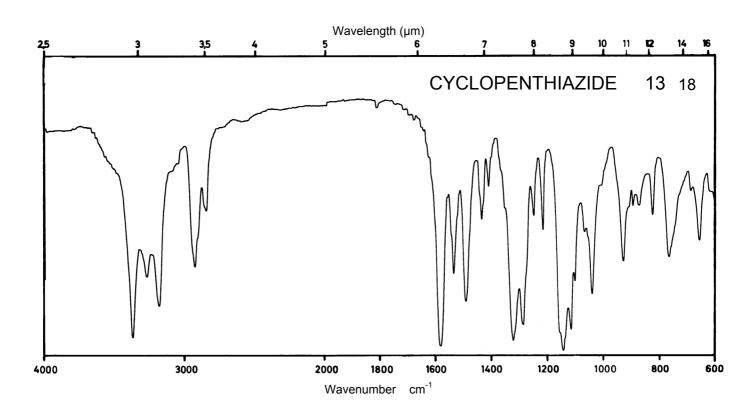



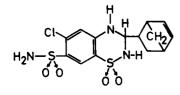



 M_{r}

357.8

Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm 270 nm		314 nm 272 nm	320 nm 274 nm
E ^{1%} 1cm	78 592		82 549	69 454
ε	2790 21180		2930 19640	2470 16240

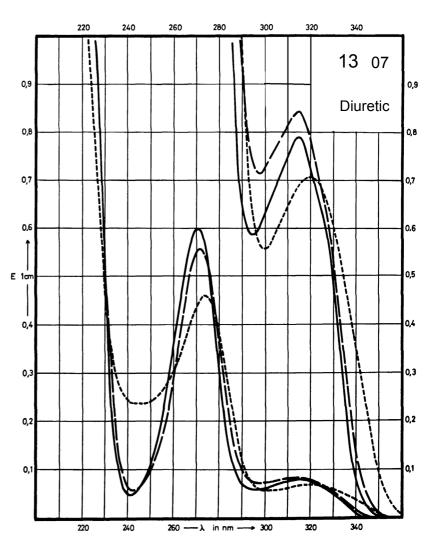

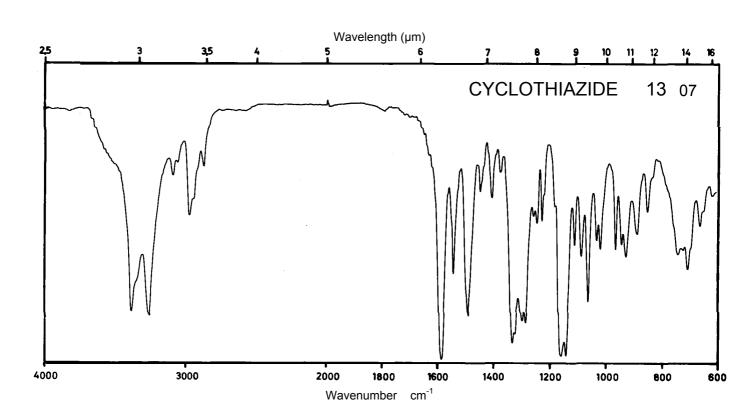



Name CLOPAMIDE

 M_r 345.9

Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	286 nm 277 nm		240 nm	
E ^{1%} 1cm	36 48		346	
ε	1250 1660		11970	

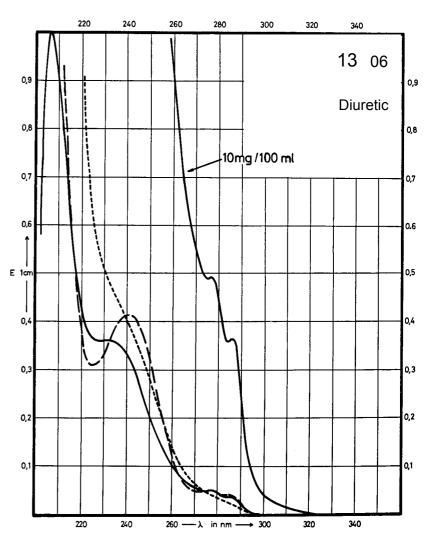

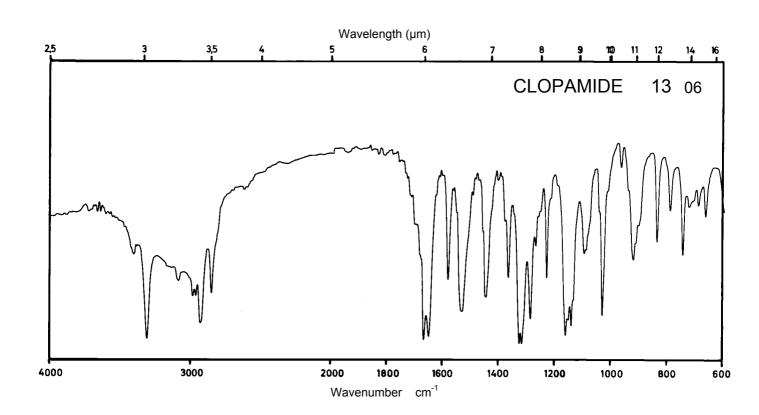

Name CHLOROTHIAZIDE


 M_{r}


295.7

Concentration


1.1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm		278 nm	291 nm
E ^{1%} 1cm	360		408	445
ε	10650		12060	13160

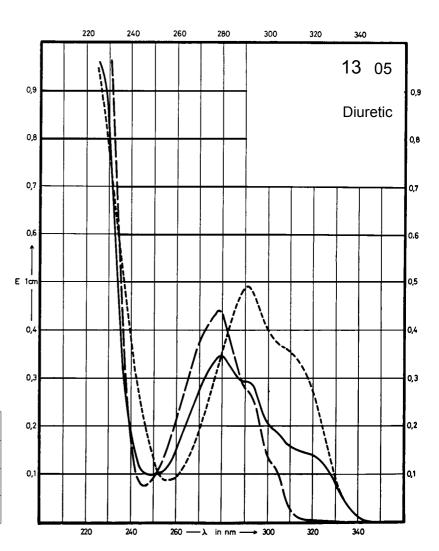

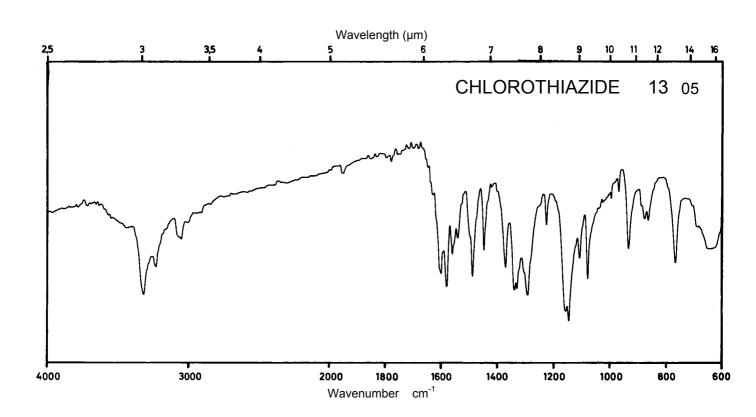



Name CHLORAZANIL HYDROCHLORIDE

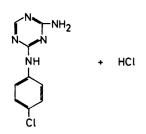

 M_r 258.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm		266 nm	264 nm
E ^{1%} 1cm	616		430	550
3	15910		11110	14210

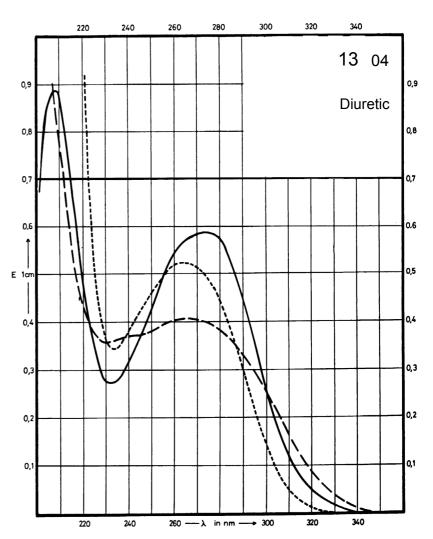

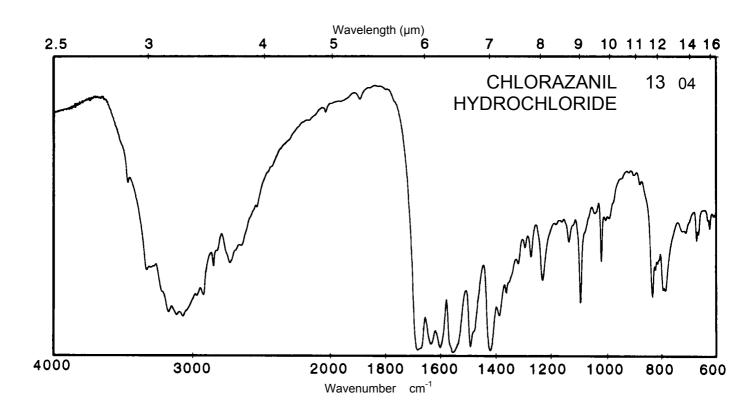


Name CHLORTALIDONE


OH CI S-NH₂

 M_{r}


338.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm 274 nm		283 nm 275 nm	
E ^{1%} 1cm	43 54		52 65	
ε	1460 1830		1760 2200	

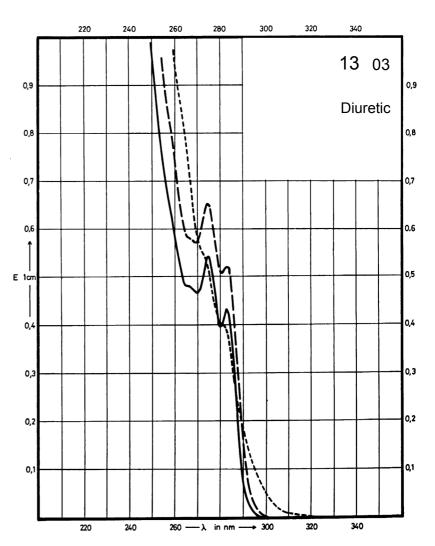

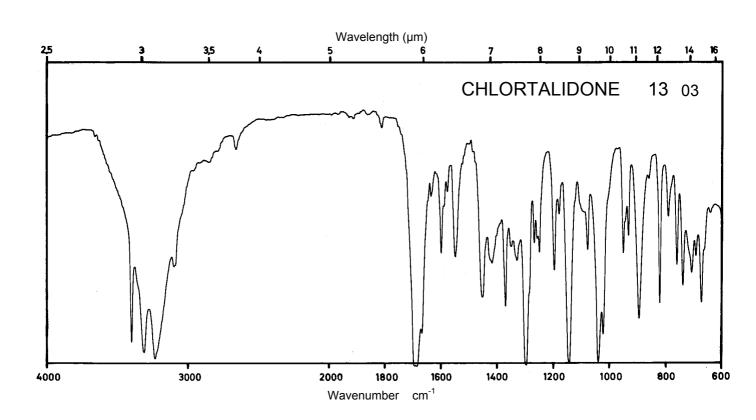




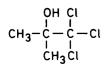


 M_{r}

177.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

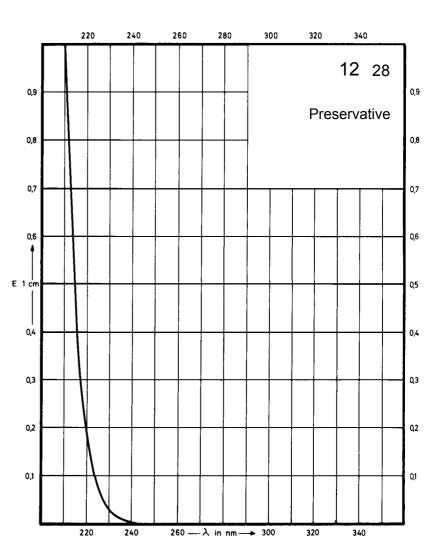


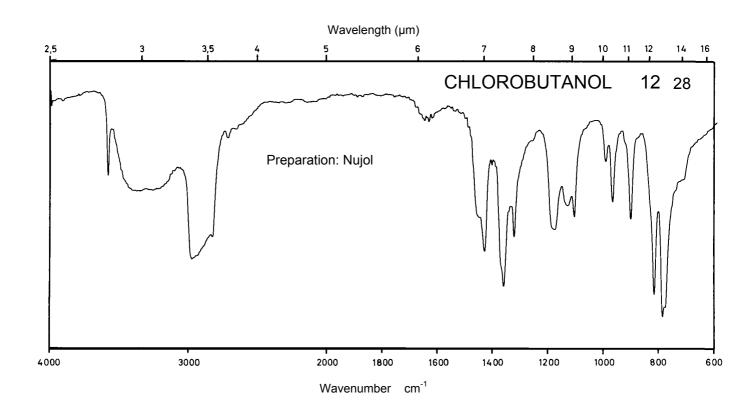


Name CHLORHEXIDINE DIHYDROCHLORIDE

 M_r 578.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm			240 nm
E ^{1%} 1cm	695			441
3	40200			25500

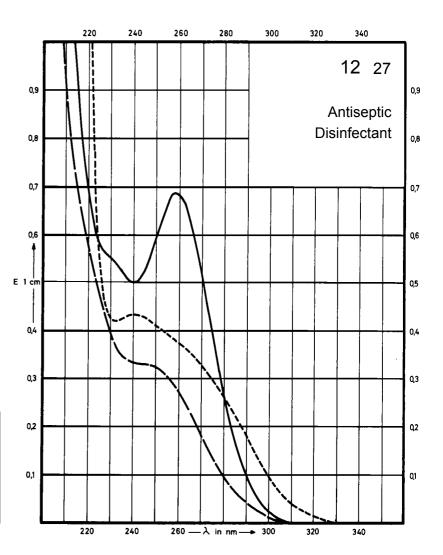


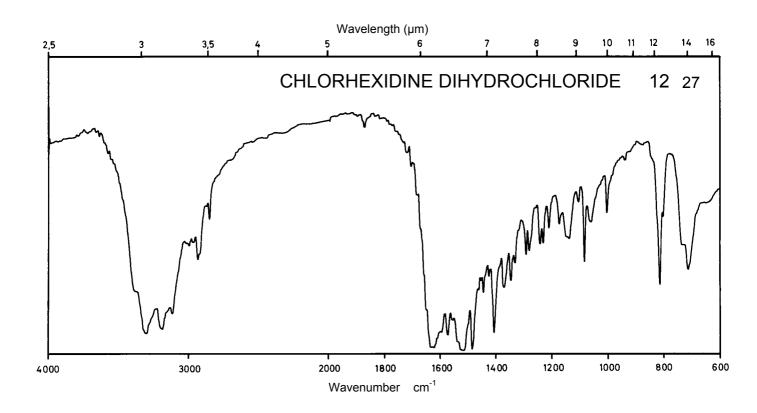


Name CHLORHEXIDINE DIACETATE

 M_r 625.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm	255 nm 232 nm		238 nm
E ^{1%} 1cm	615	446 454		366
3	38480	27930 28410		22930

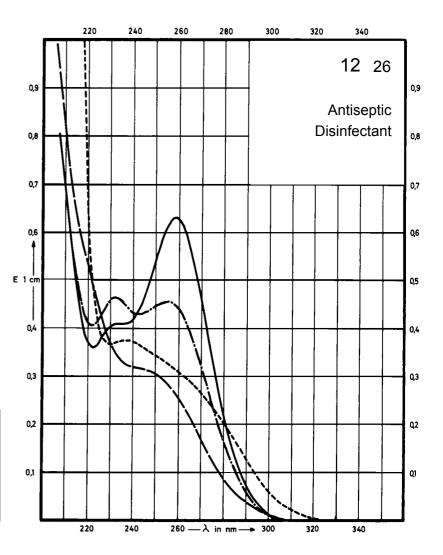

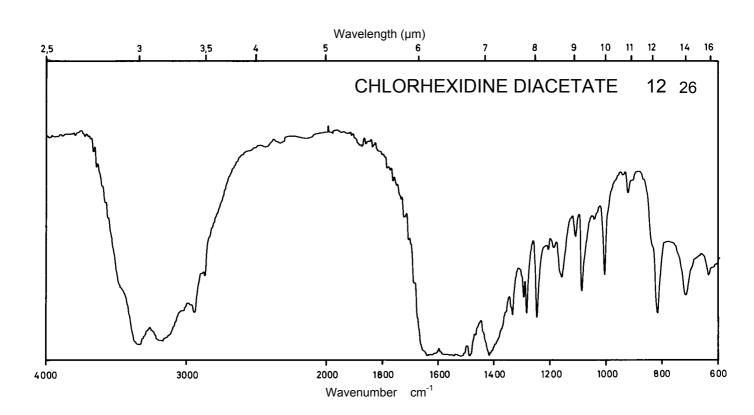



Name CETYLPYRIDINIUM CHLORIDE

 M_{r}


358.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm		260 nm	260 nm
E ^{1%} 1cm	126		121	122
ε	4510		4330	4370

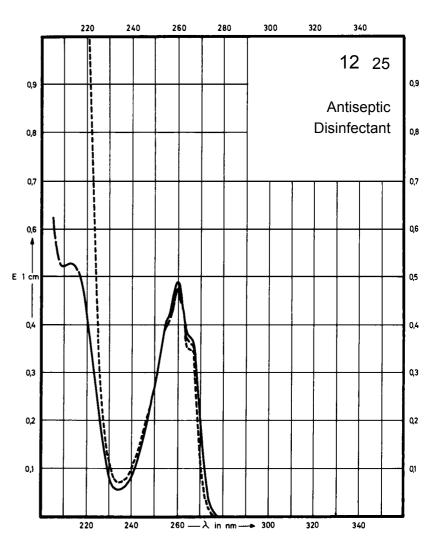

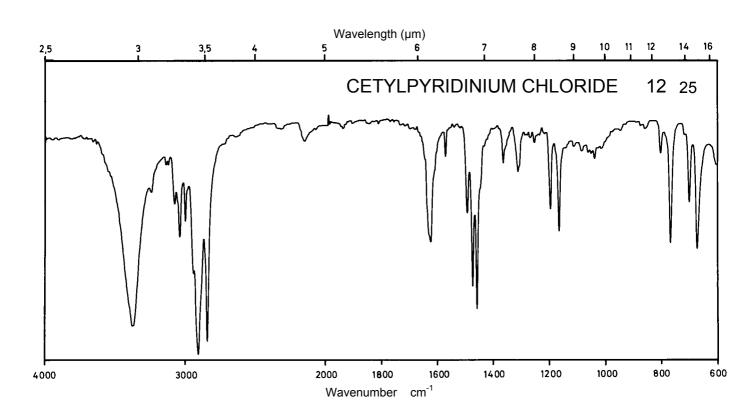



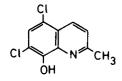



228.1

Concentration

0.25 mg / 100 ml 0.5 mg / 100 ml 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm 251 nm		330 nm 263 nm	345 nm 263 nm
E ^{1%} 1cm	146 1900		132 2100	186 1410
3	3340 43340		3010 47900	4240 32160

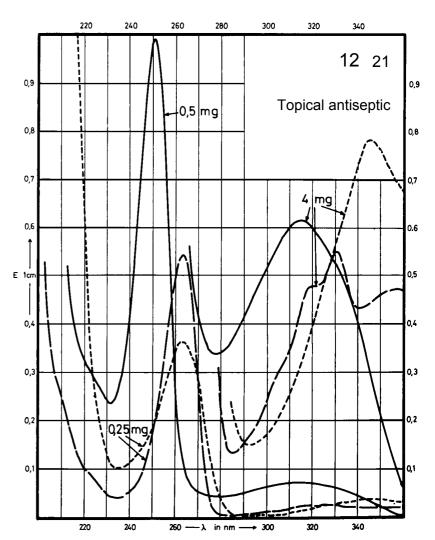

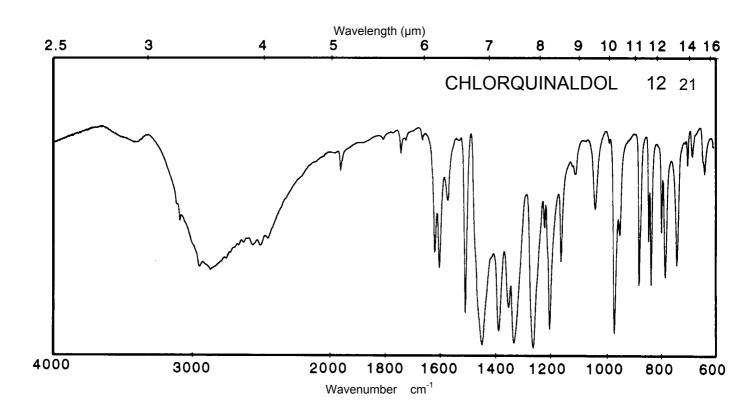



Name 4-CRESOL

OH CH₃

 M_r 108.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	279 nm 223 nm		277 nm 220 nm	295 nm 236 nm
E ^{1%} 1cm	186 610		162 550	245 800
3	2010 6590		1750 5950	2650 8650

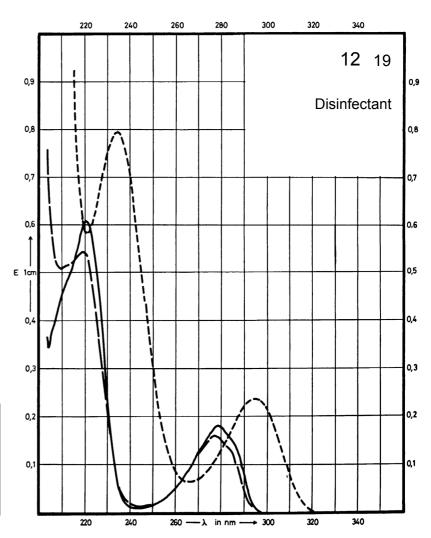

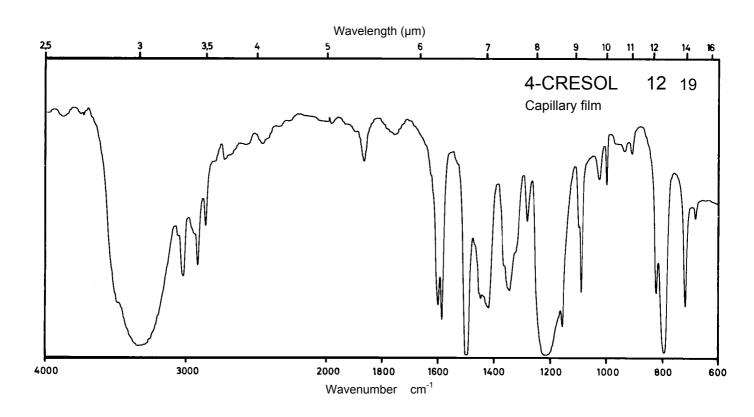


Name 2-CRESOL


OH CH3


 M_r

108.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm 214 nm		271 nm	287 nm 237 nm
E ^{1%} 1cm	190 600		155	296 840
ε	2050 6490		1680	3200 9080

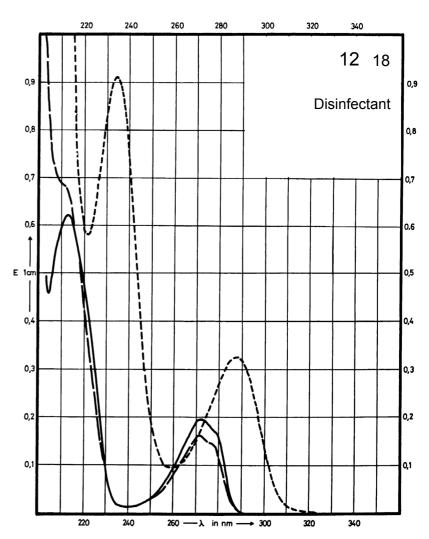

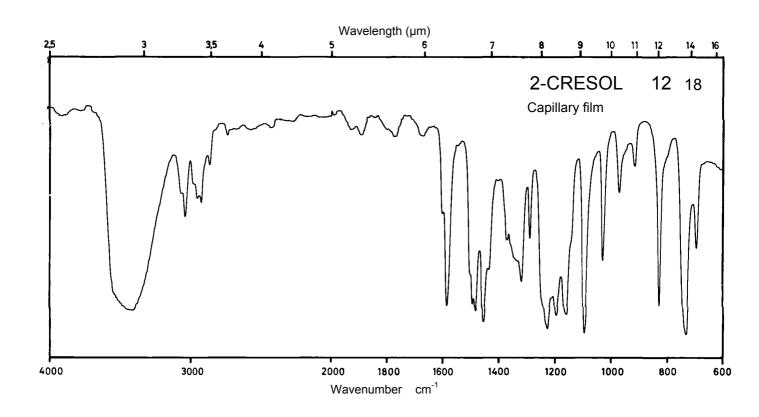


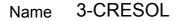





 M_r 108.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm		271 nm	288 nm 237 nm
E ^{1%} 1cm	165		139	251 825
3	1785		1505	2710 8920

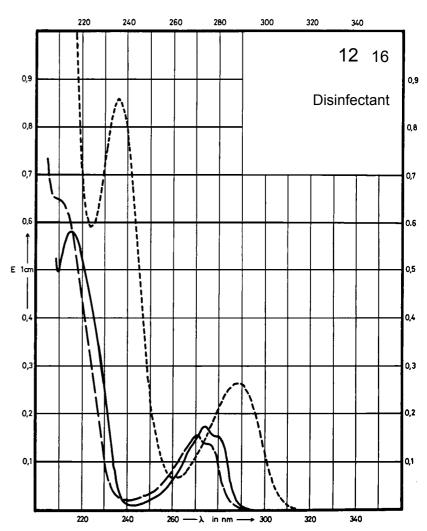

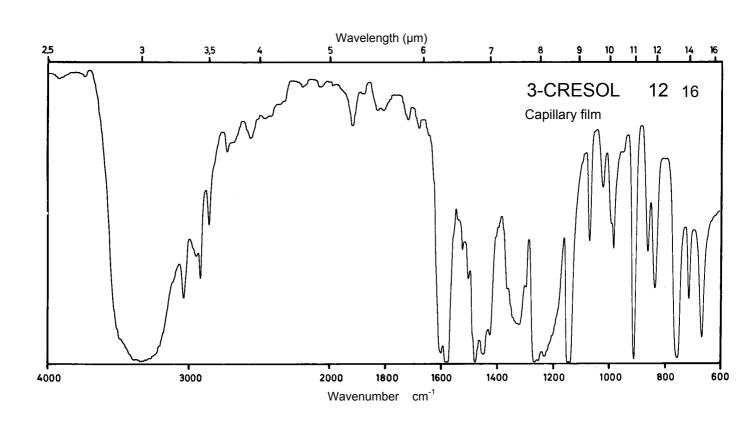


Name CHLOROXYLENOL

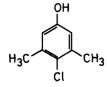


 M_{r}

156.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm 219 nm		277 nm 218 nm	296 nm 242 nm
E ^{1%} 1cm	97 433		84 460	156 650
ε	1520 6780		1320 7200	2440 10180

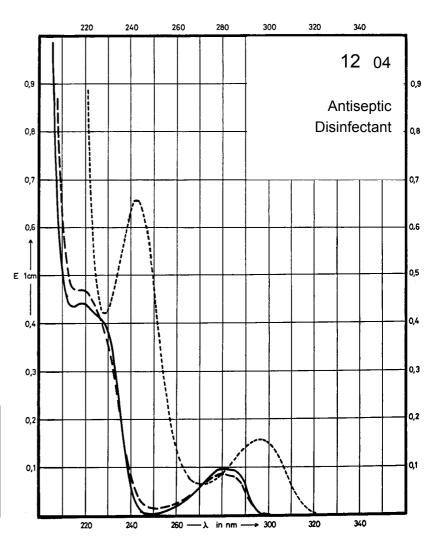

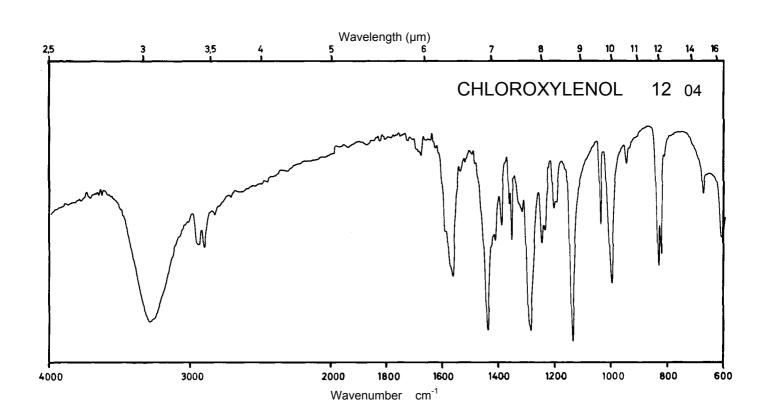

Name CHLOROCRESOL

OH CI


 M_{r}


142.6

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 228 nm		278 nm 225 nm	296 nm 241 nm
E ^{1%} 1cm	127 534		107 497	169 732
ε	1810 7610		1530 7090	2410 10440

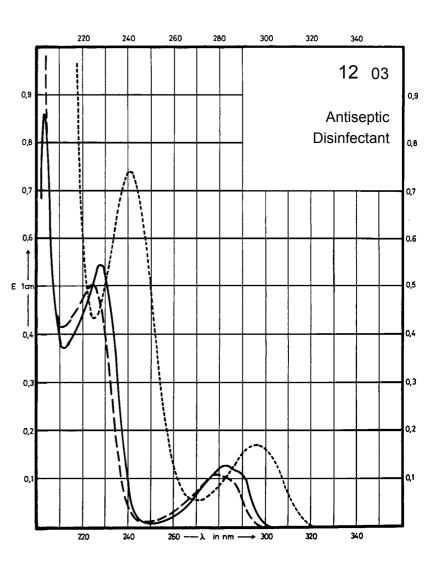

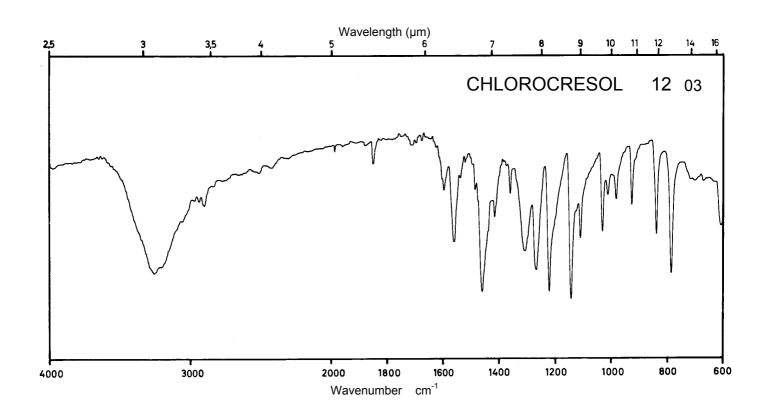





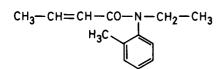


203.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

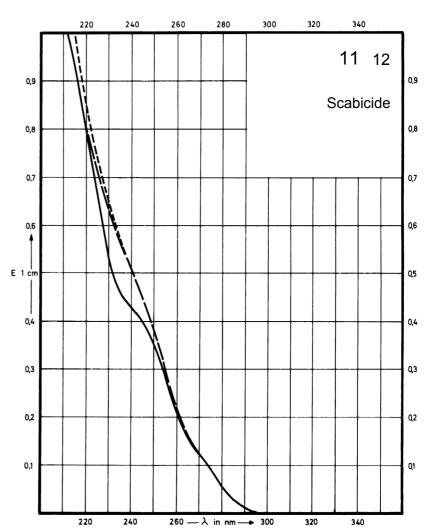

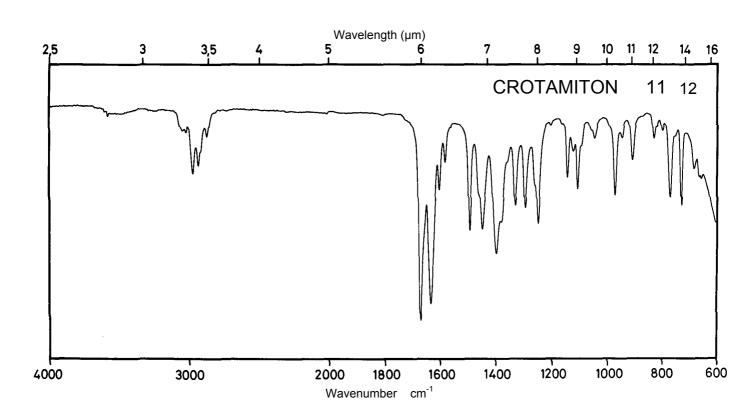





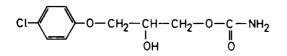


 M_{r} 245.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 281 nm		279 nm	277 nm
E ^{1%} 1cm	53 66		57	57
3	1300 1620		1400	1400

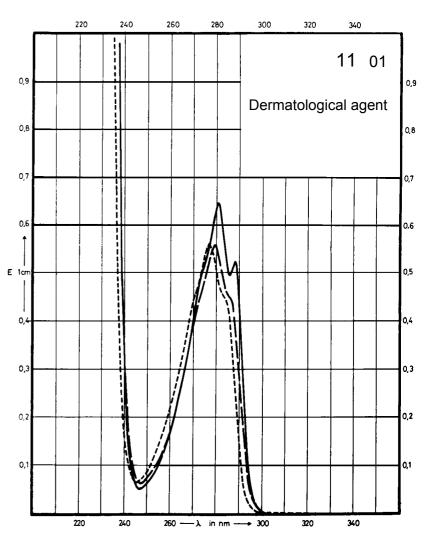

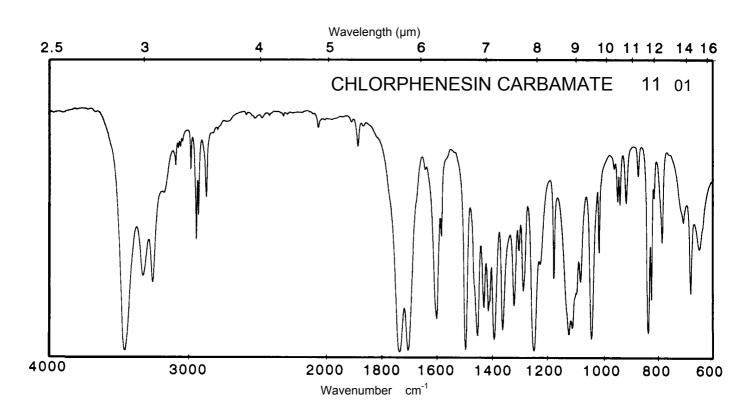


Name CLOPREDNOL

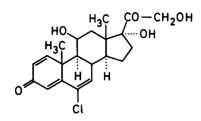


 M_{r}

392.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	298 nm 255 nm 227 nm	306 nm 260 nm 227 nm		
E ^{1%} 1cm	265 265 284	251 269 319		
ε	10400 10400 11150	9880 10600 12500		

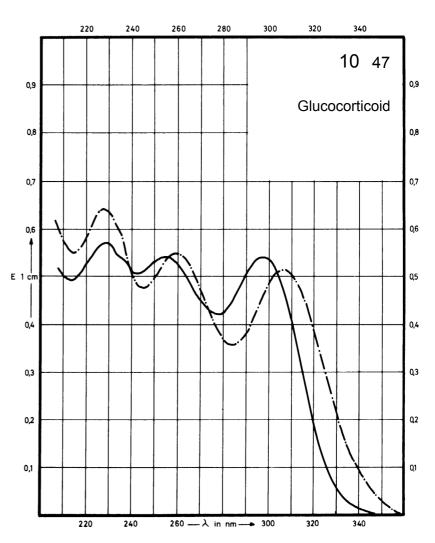

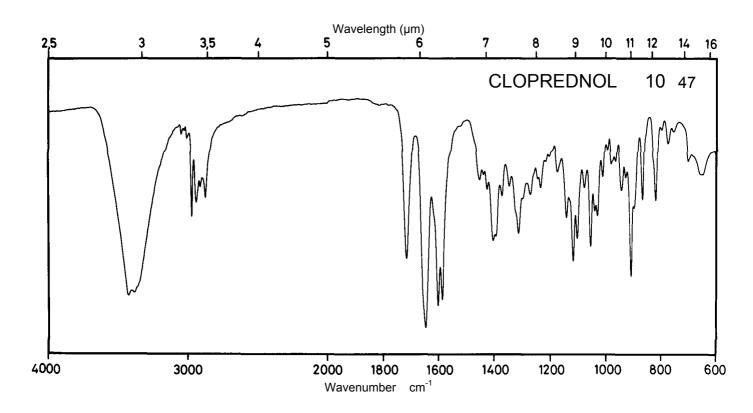



Name CLOCORTOLONE 21-PIVALATE

 M_r


495.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	237 nm			
E ^{1%} 1cm	328			
3	16250			

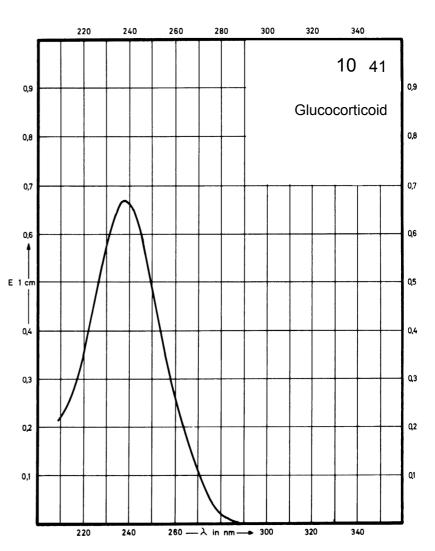

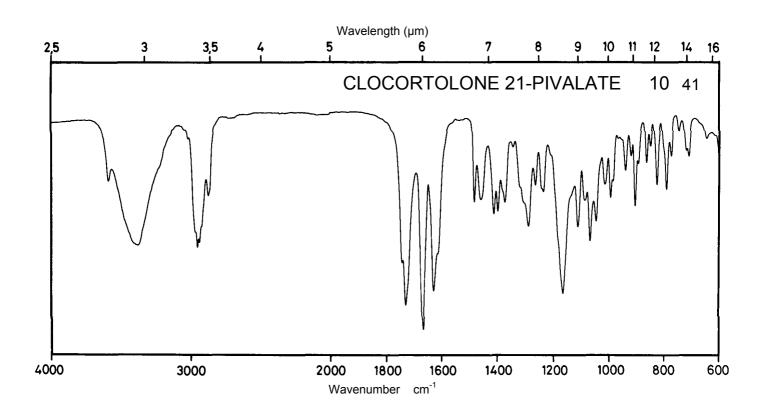



Name CLOBETASOL 17α-PROPIONATE

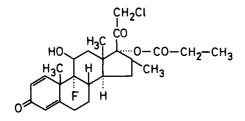


 M_{r} 467.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm	242 nm		
E ^{1%} 1cm	339	332		
ε	15800	15500		

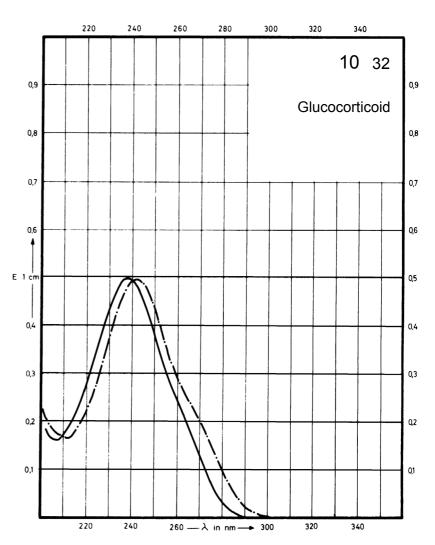

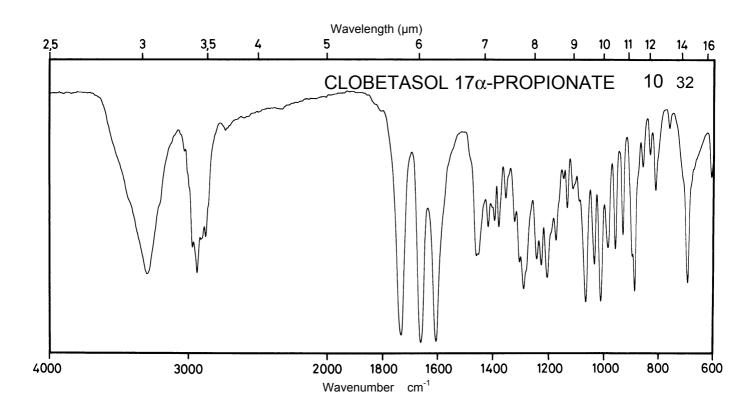





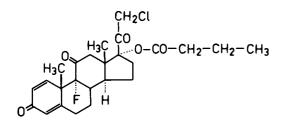


479.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	236 nm			
E ^{1%} 1cm	330			
ε	15800			

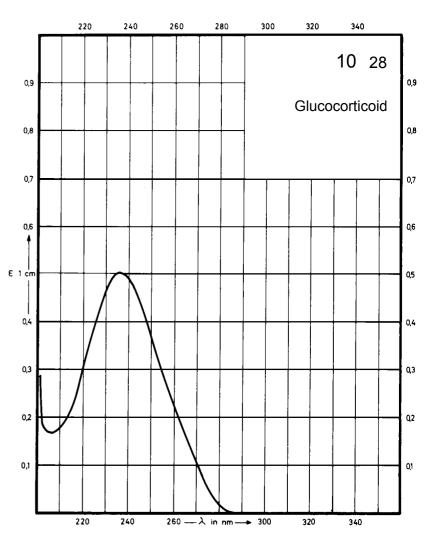

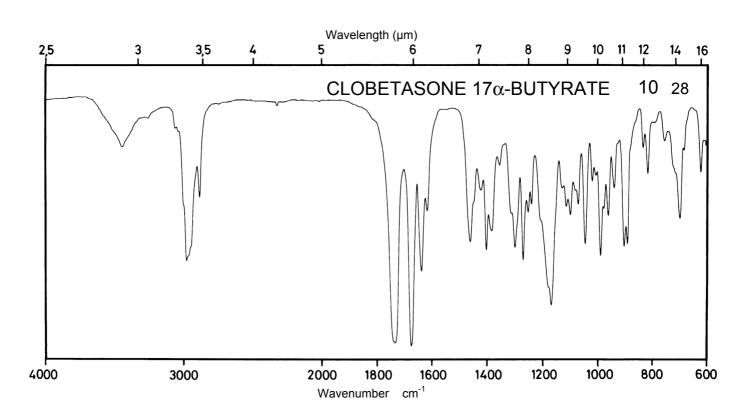





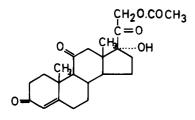



402.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	237 nm			
E ^{1%} 1cm	395			
3	15900			

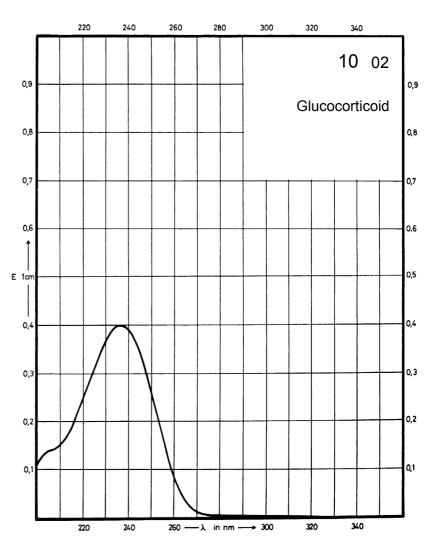

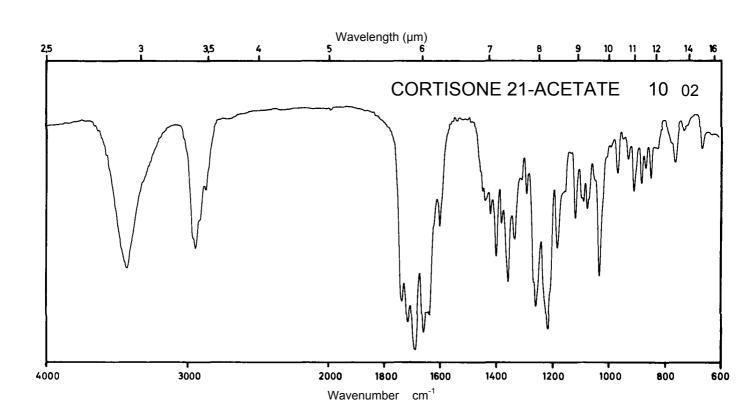

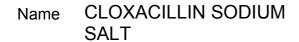


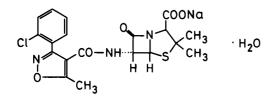


 M_r 475.9


Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm 266 nm		Decom- position observed	Decom- position observed
E ^{1%} 1cm	8.5 9.9			
ε	405 470			

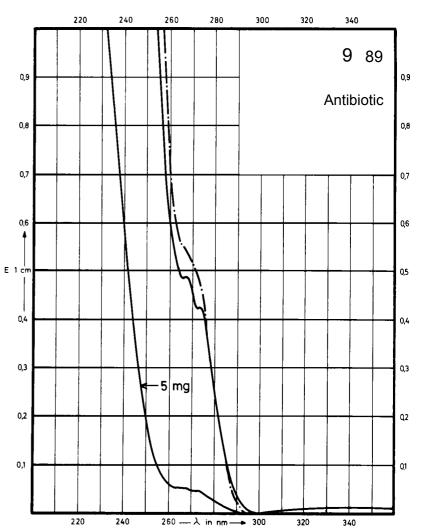

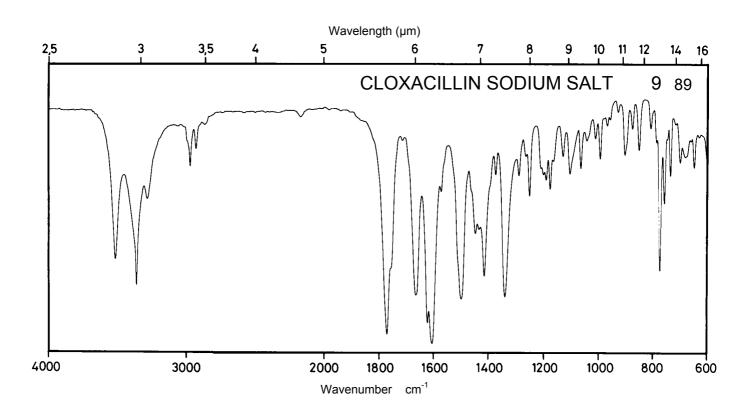




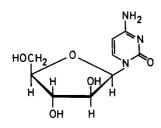


 M_r 243.2

Concentration 1 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		280 nm	274 nm
E ^{1%} 1cm	372		560	416
3	9050		13600	10100

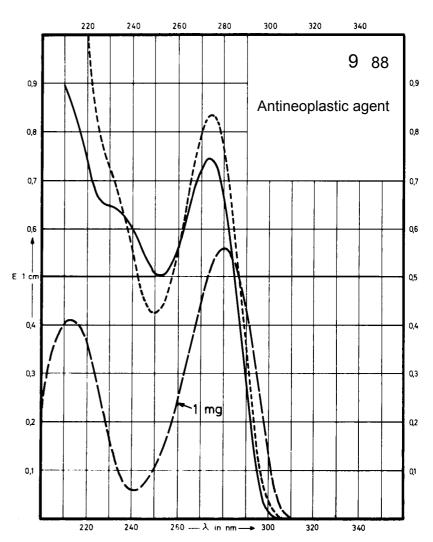

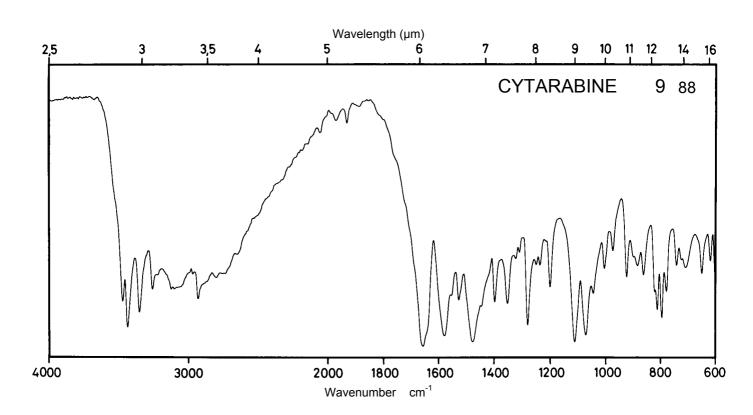







102.1

Concentration


2.2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		226 nm	218 nm	
E ^{1%} 1cm		380	345	
ε		3900	3500	

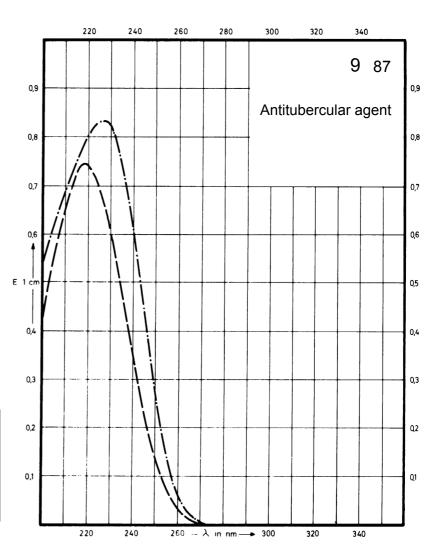

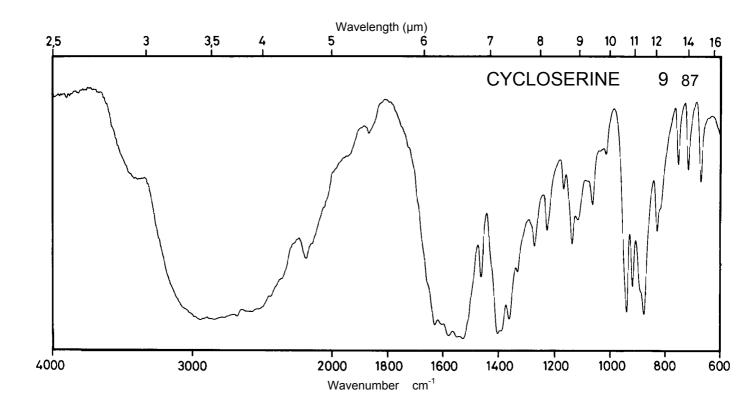





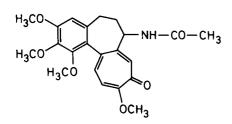


399.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	350 nm 242 nm		354 nm 247 nm	353 nm 245 nm
E ^{1%} 1cm	427 775		411 870	430 851
ε	17100 31000		16400 34700	17200 34000

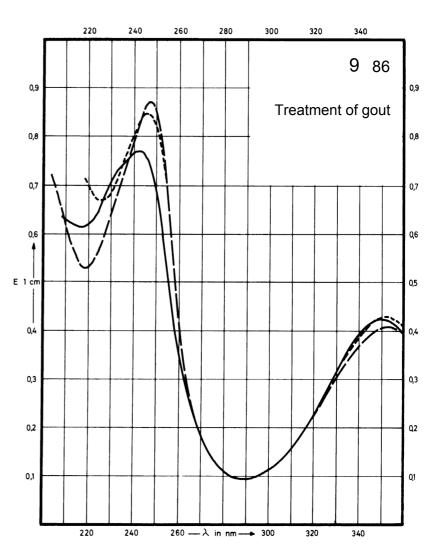

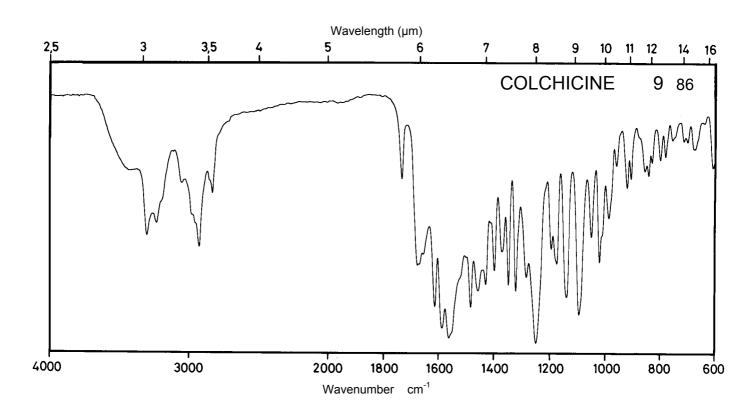




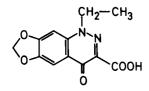


 M_r 262.2

Concentration 0.5 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	352 nm 257 nm		354 nm 263 nm	351 nm 268 nm 244 nm
E ^{1%} 1cm	488 950		446 1100	547 787 808
ε	12800 24900		11700 28800	14300 20600 21200

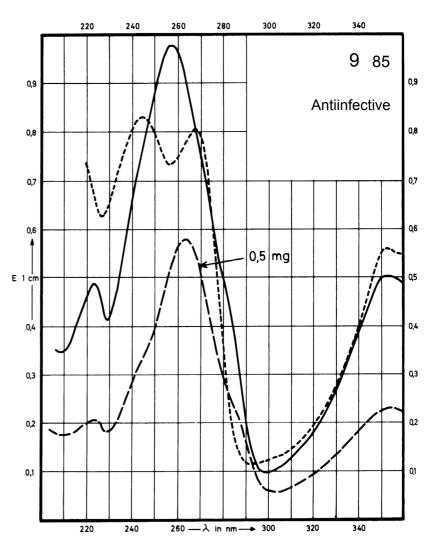

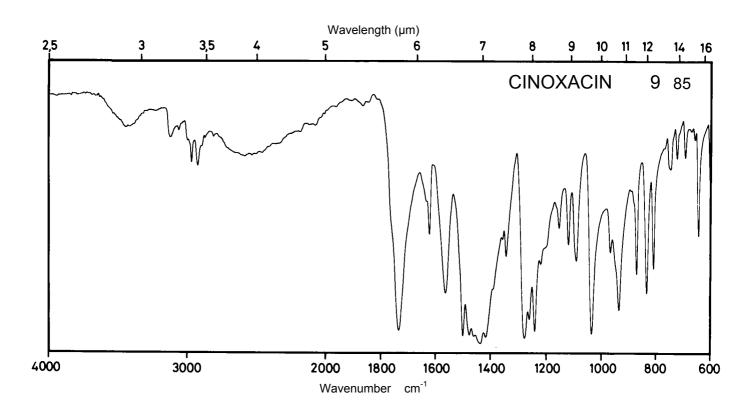

Name CLOXIQUINE

OH N


 M_{r}


179.6

Concentration

0.4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	329 nm 245 nm		381 nm 256 nm	368 nm 255 nm
E ^{1%} 1cm	194 2090		121 2380	205 1600
3	3480 37540		2170 42740	3680 28740

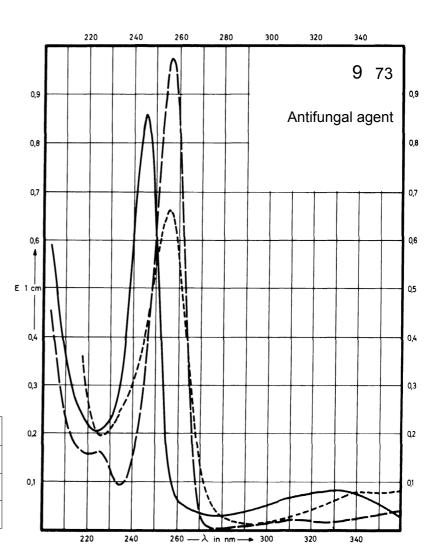

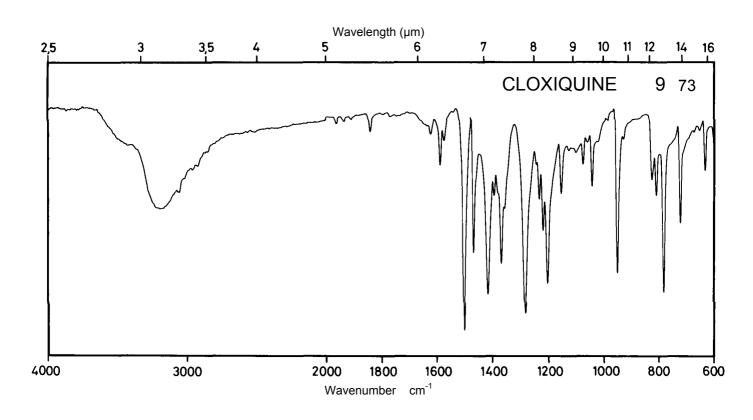







 M_r 351.1


Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 253 nm		262 nm	263 nm
E ^{1%} 1cm	97 1190		830	1055
3	3390 41800		29140	37040

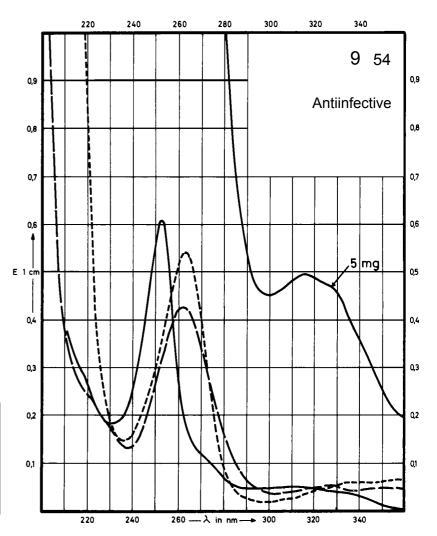

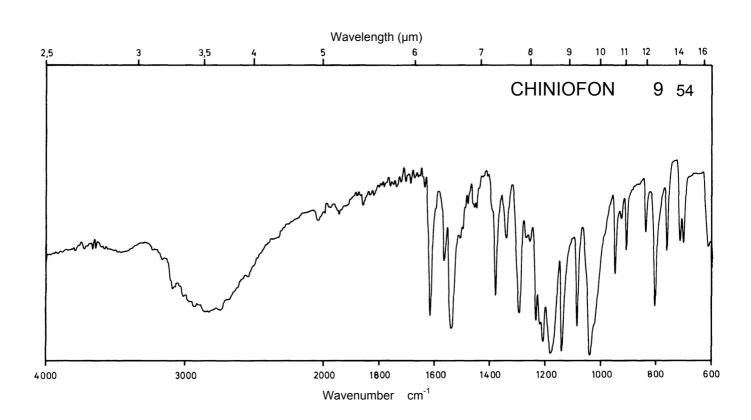





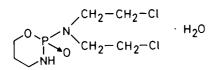

279.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

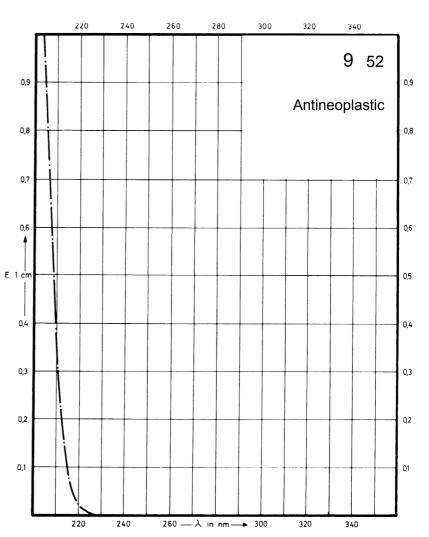

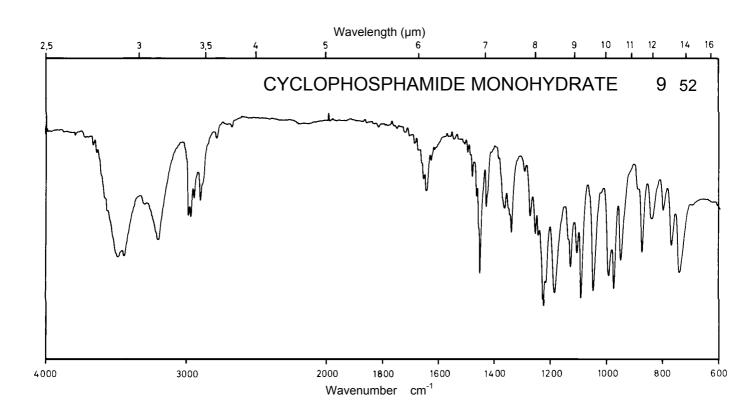





$$\begin{array}{c} \text{CI-CH}_2\text{-CH}_2\\ \text{CI-CH}_2\text{-CH}_2\end{array} \text{N-CH}_2\text{-CH}_2\text{-CH}_2\text{-COOH} \\ \end{array}$$

 M_{r} 304.2

Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 303 nm 258 nm 303 nm 258 nm 258 nm absorption E^{1%} 1cm 67 634 32 2050 19280 1800 15900 ε 970

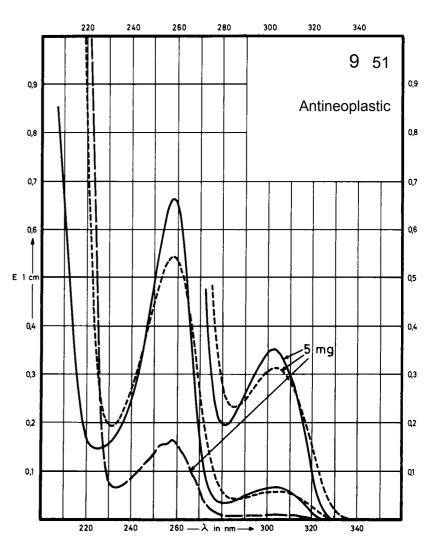

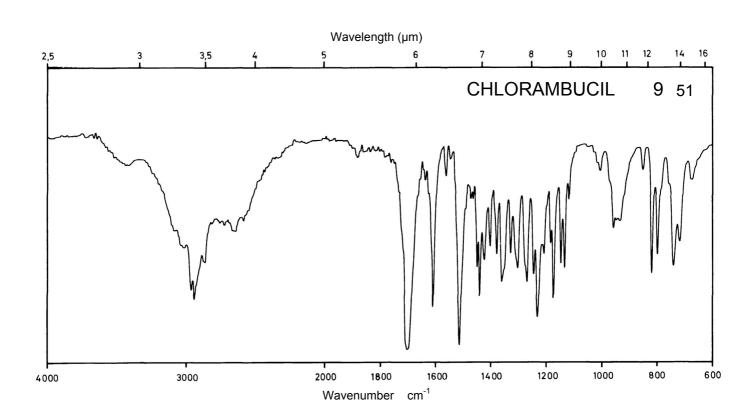


Name CLOTRIMAZOLE




 M_{r}

344.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm		262 nm	
E ^{1%} 1cm	22		26	
ε	760		900	

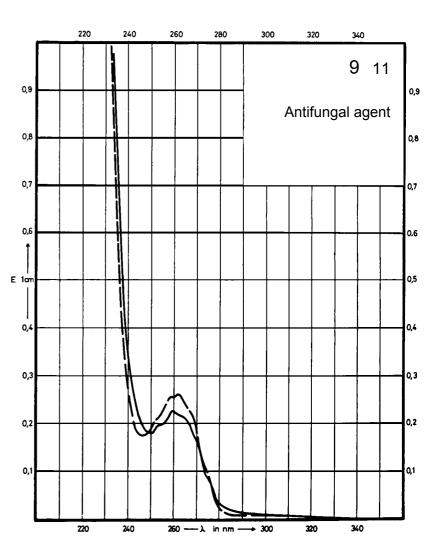

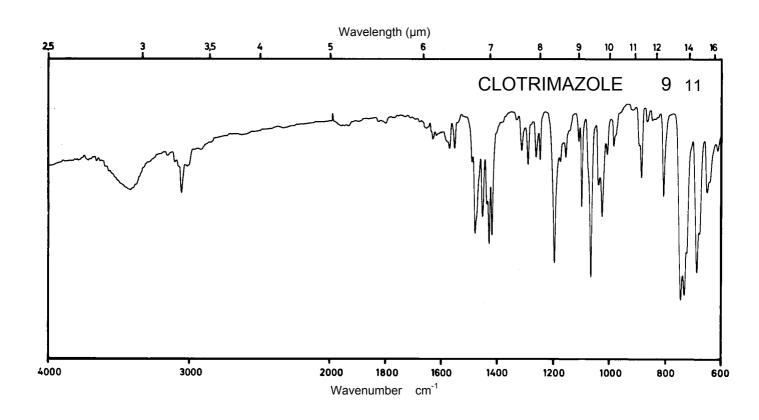







305.5


Concentration 0.5 mg / 100 ml 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	323 nm 254 nm		333 nm 264 nm	350 nm 265 nm
E ^{1%} 1cm	109 1133		91 952	156 1154
ε	3330 34610		2780 29080	4770 35250

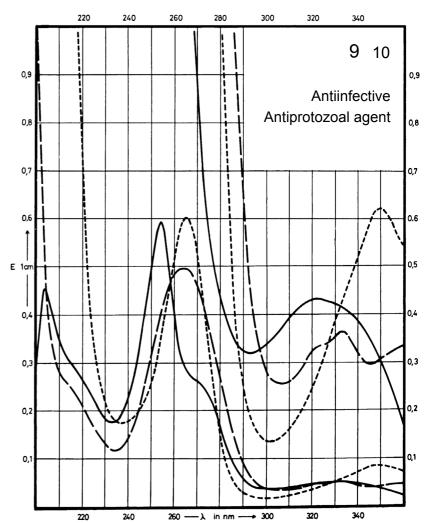

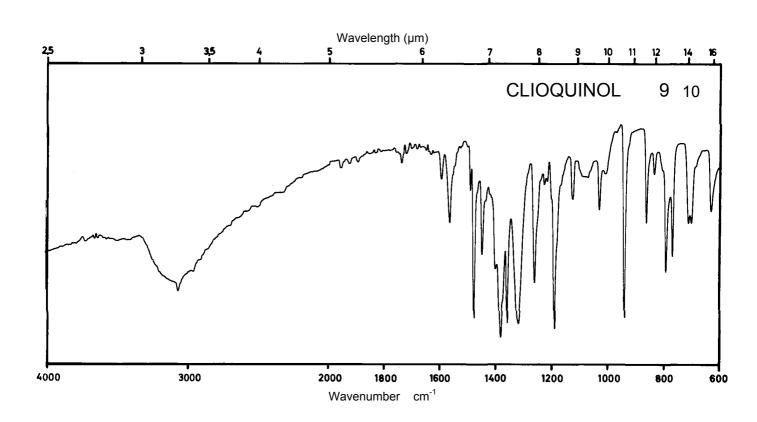


Name CICLOPIROX OLAMINE

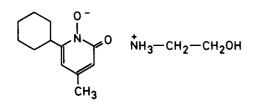


 M_{r}

268.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	303 nm		292 nm	312 nm
E 1% 1cm	264		272	257
ε	7090		7300	6900

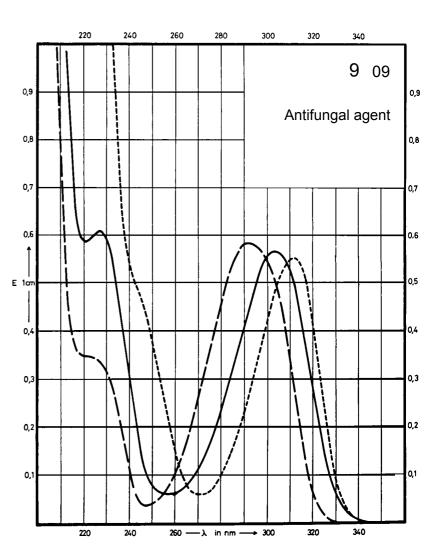

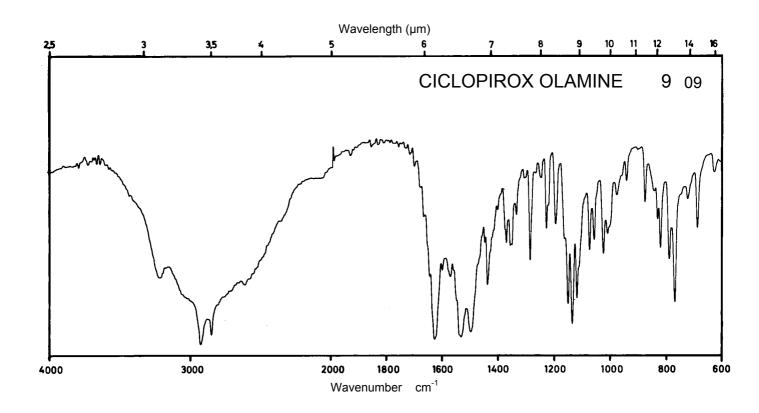



Name CHLOROQUINE PHOSPHATE

 M_{r}


515.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	330 nm 254 nm		343 nm 257 nm	330 nm 254 nm
E ^{1%} _{1cm}	268 337		378 325	220 344
ε	13830 17390		19500 16770	11350 17750

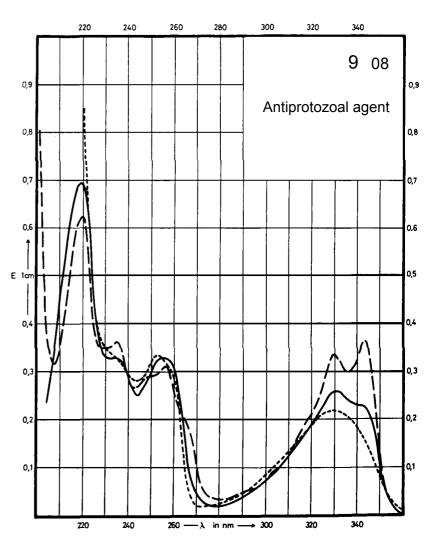

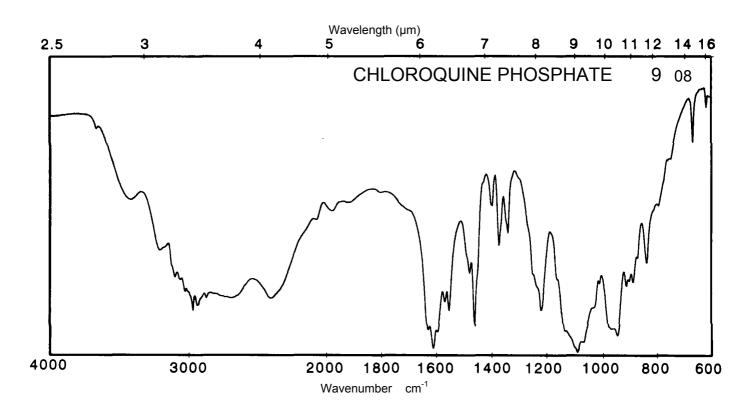


Name CLOFEDANOL HYDROCHLORIDE

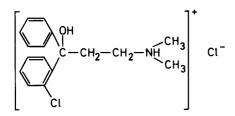



 M_{r}

326.3


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 259 nm		264 nm 258 nm	
E ^{1%} 1cm	10.0 10.8		10.4 11.2	
3	326 352		340 365	

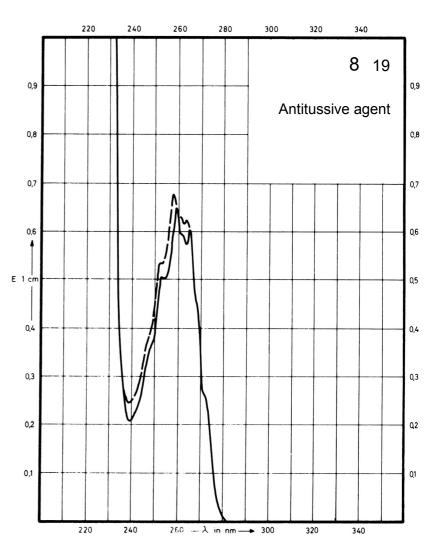


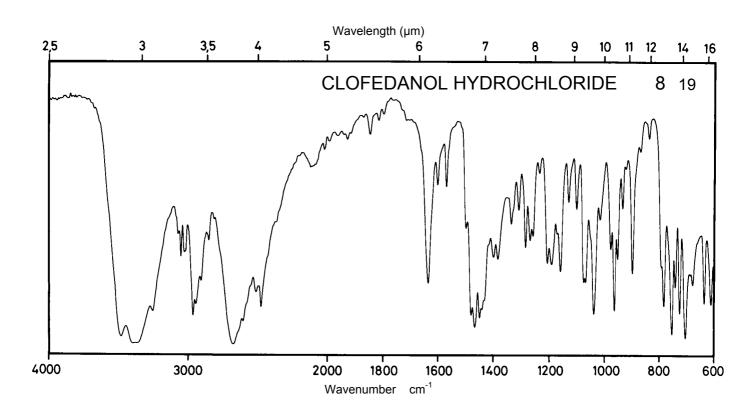

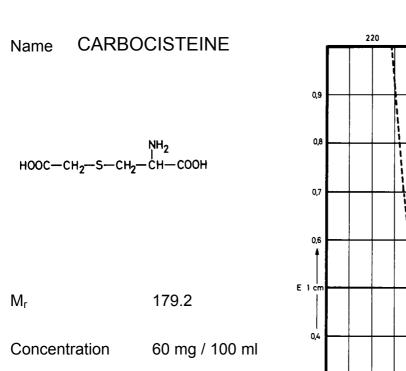


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
3				

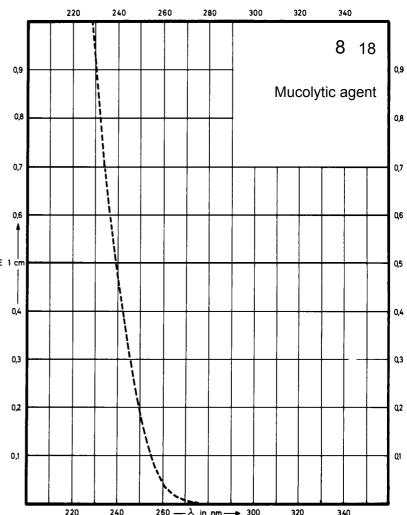

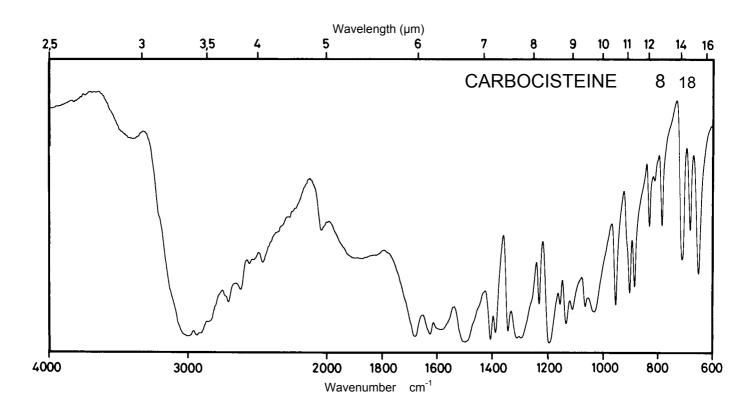


Name CLOBUTINOL HYDROCHLORIDE

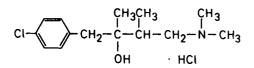


 M_{r}


292.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 268 nm	276 nm 268 nm	276 nm 268 nm	
E ^{1%} 1cm	7.5 9.8	7.0 9.7	7.0 9.7	
ε	219 286	205 283	205 283	

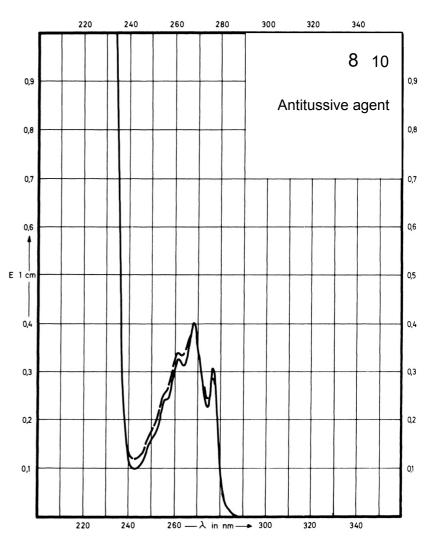

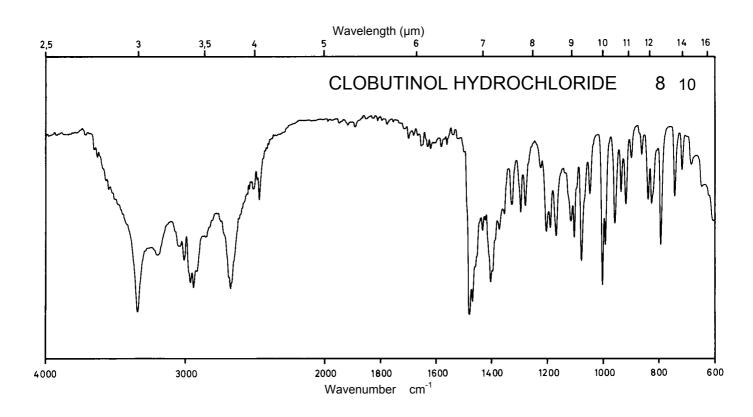





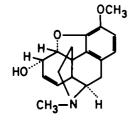


 M_{r} 299.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	286 nm			
E ^{1%} 1cm	52.8			
3	1580			

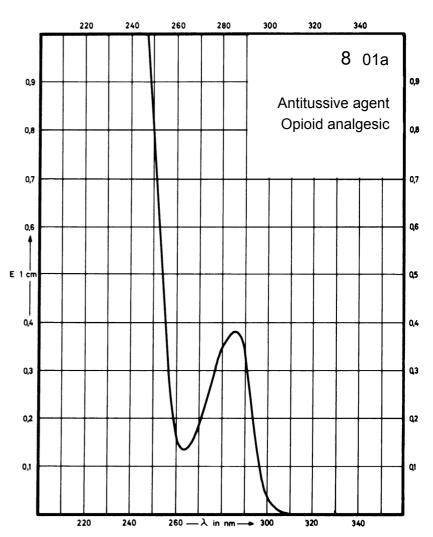

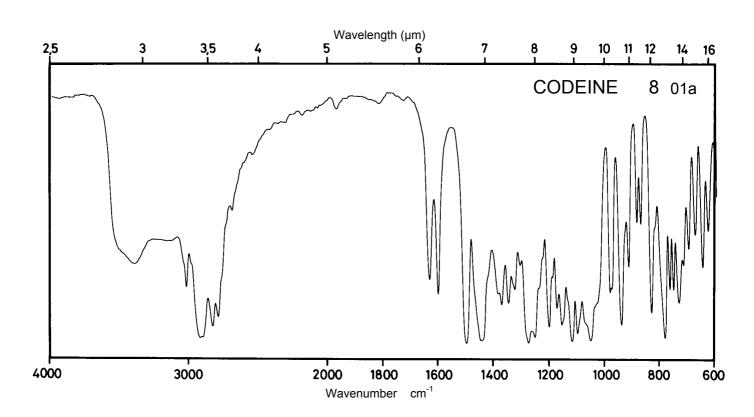



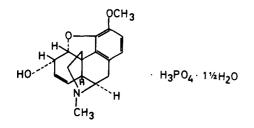



424.4


Concentration


10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm		284 nm	284 nm
E ^{1%} _{1cm}	39		38	35
ε	1660		1610	1490

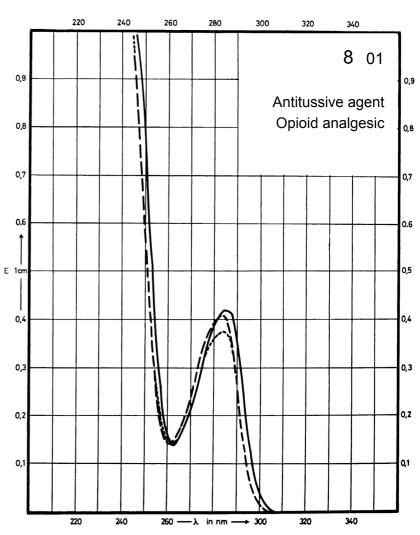

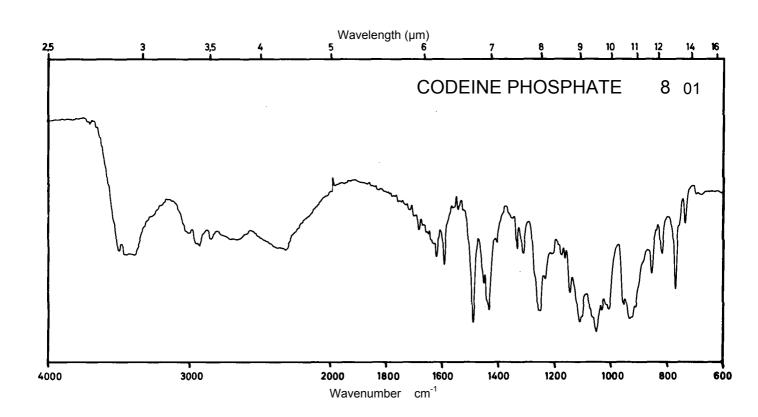

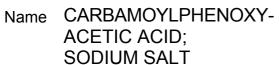

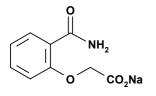



 M_r

Concentration 1 mg / 100 ml

217.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 232 nm	291 nm 235 nm	289 nm 233 nm	292 nm
E ^{1%} 1cm	144 358	137 344	121 347	140
٤	3130	2970	2640	3040

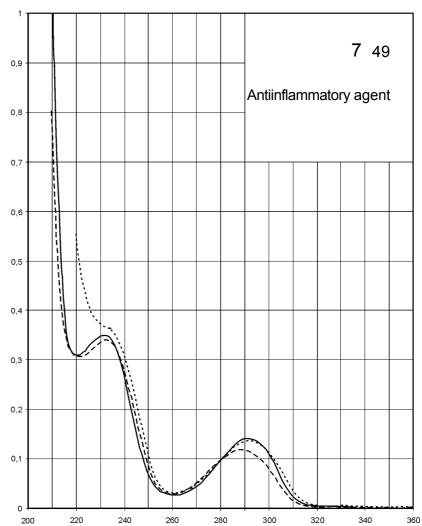

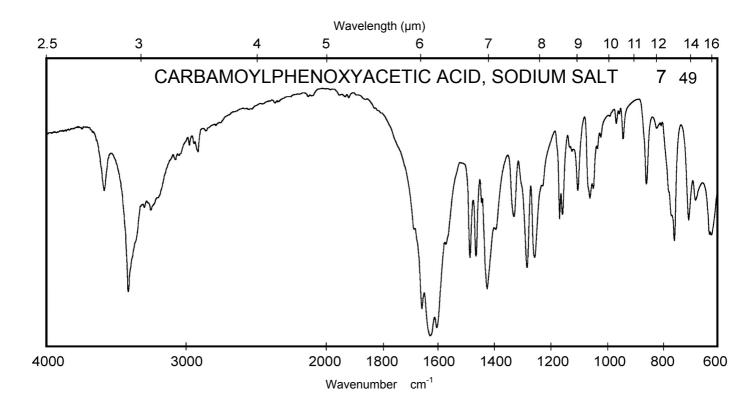




C14H21CIN2O2 · C19H20N2O2 · 2H2O


 M_r

Concentration 1.7 mg / 100 ml


629.2

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 268 nm 224 nm 225 nm 264 nm absorption E^{1%} 1cm 373 366 372 338 23500 23000 3 23400 21300

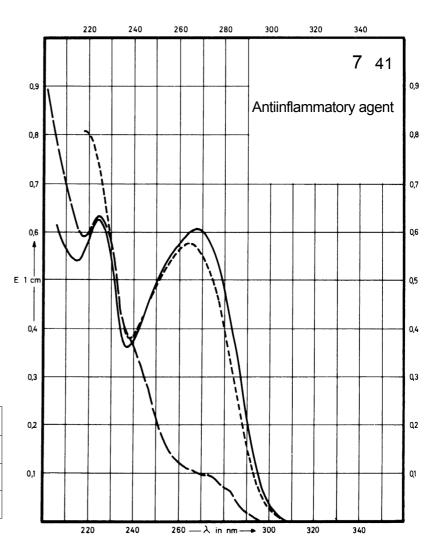

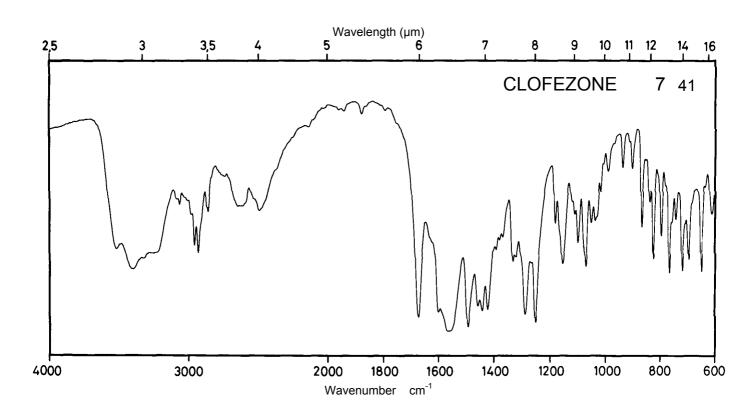



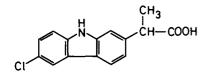



273.7

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	300 nm 262 nm 238 nm			300 nm 262 nm 238 nm
E ^{1%} 1cm	775 825 1760			702 791 1675
ε	21200 22600 48200			19200 21600 45800

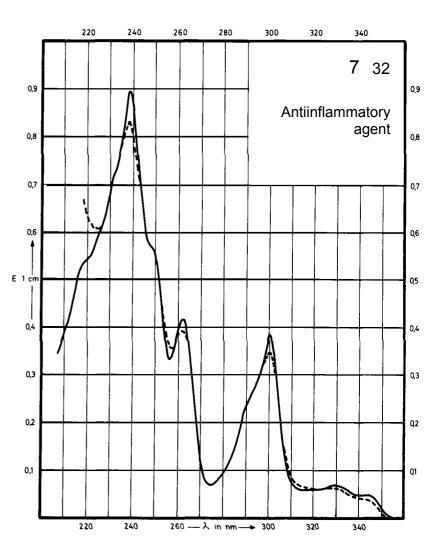

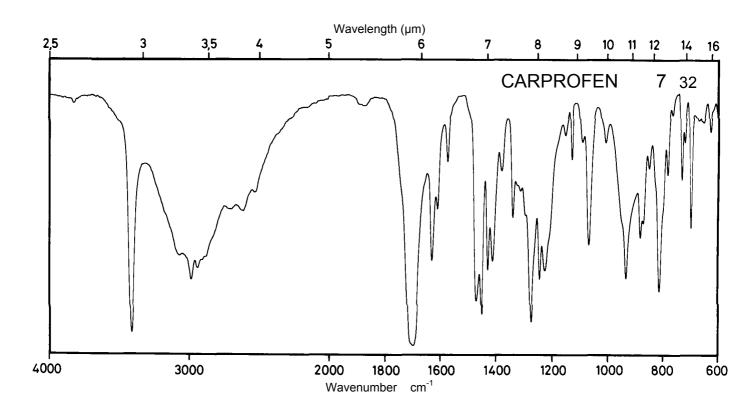




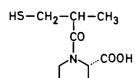


217.3


Concentration 1 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

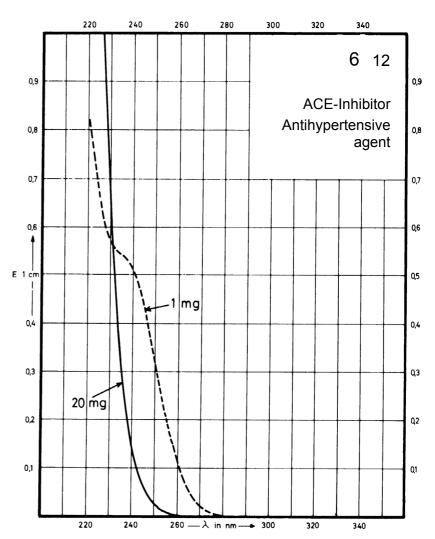

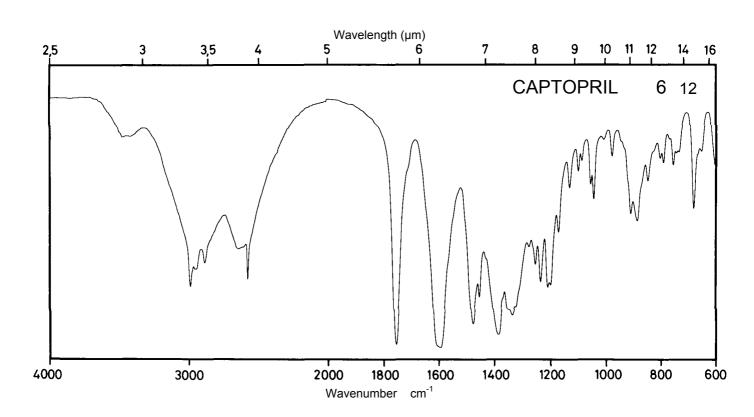


Name CLONIDINE HYDROCHLORIDE

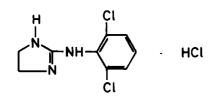


 M_{r}


266.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm 271 nm		278 nm 271 nm	
E ^{1%} 1cm	15.6 18.1		14.2 16.6	
ε	420 480		380 440	

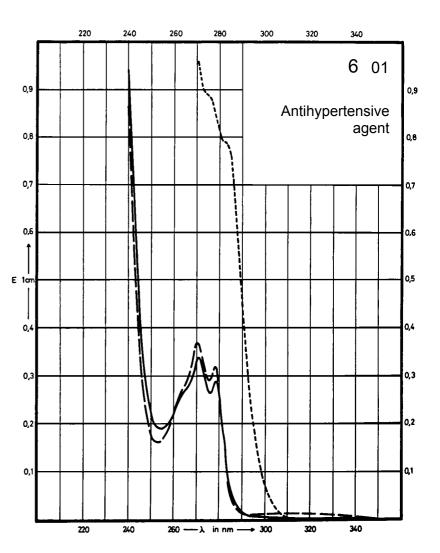

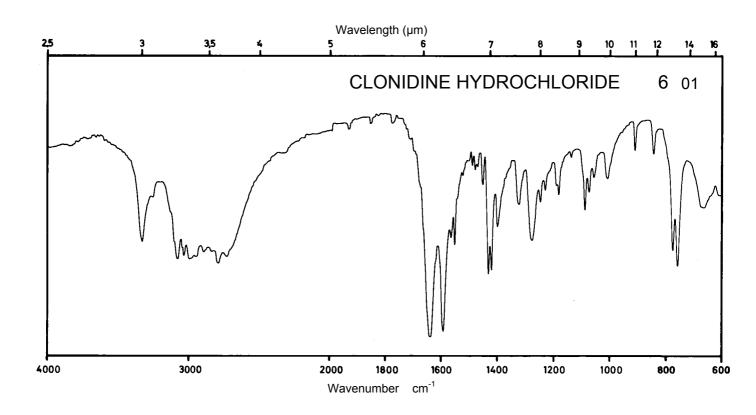




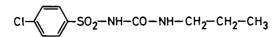


 M_r 276.7


Concentration 1 mg / 100 ml 25 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 276 nm 265 nm 231 nm 15.7 20.4 Maximum of 232 nm 228 nm absorption E^{1%} 1cm 598 472 575 430 560 15900 ε 16500 13100

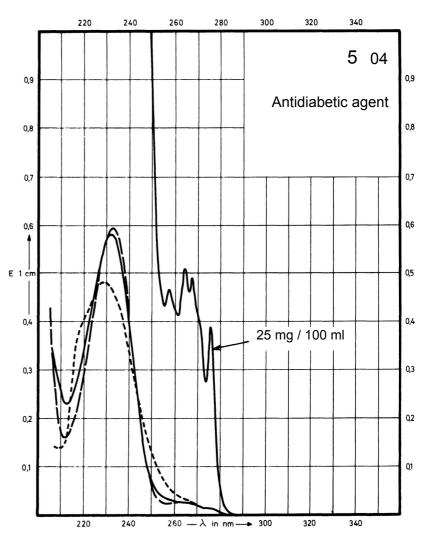

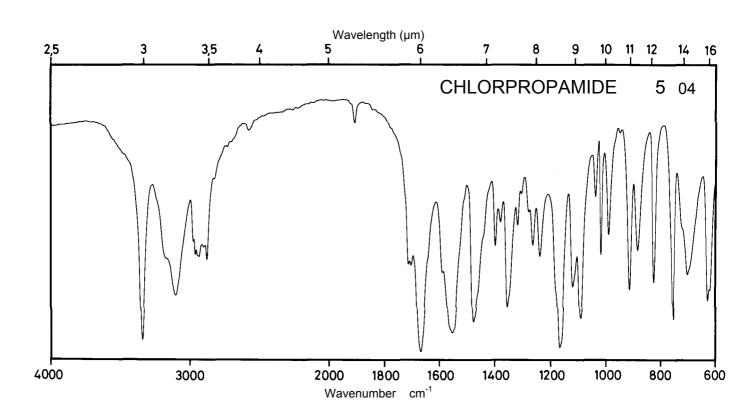




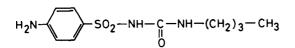


271.4


Concentration 1 mg / 100 ml 2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm		265 nm	254 nm
E ^{1%} 1cm	740		163	635
3	20080		4420	17230

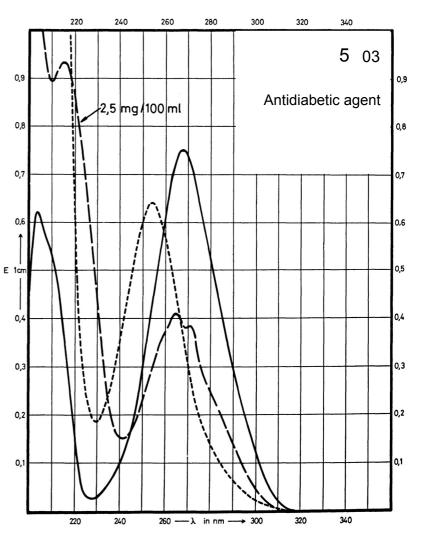

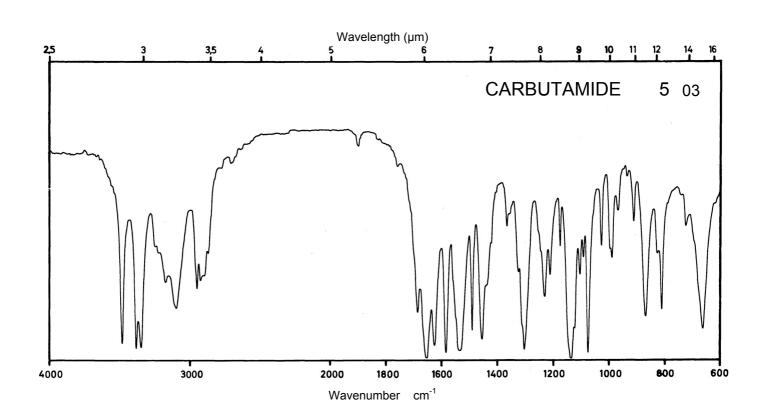



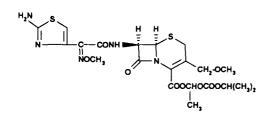



557.6

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	235 nm	232 nm	263 nm	
E ^{1%} 1cm	330	293	338	
3	18400	16300	18800	

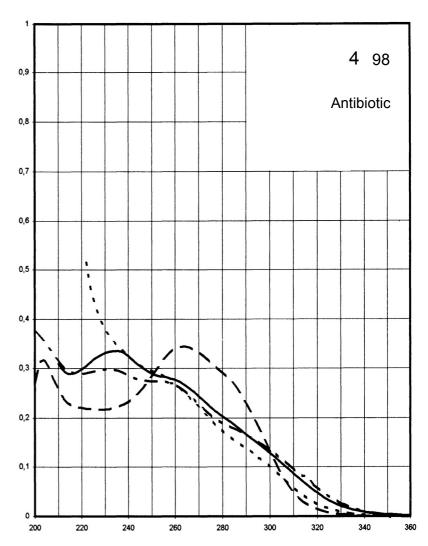

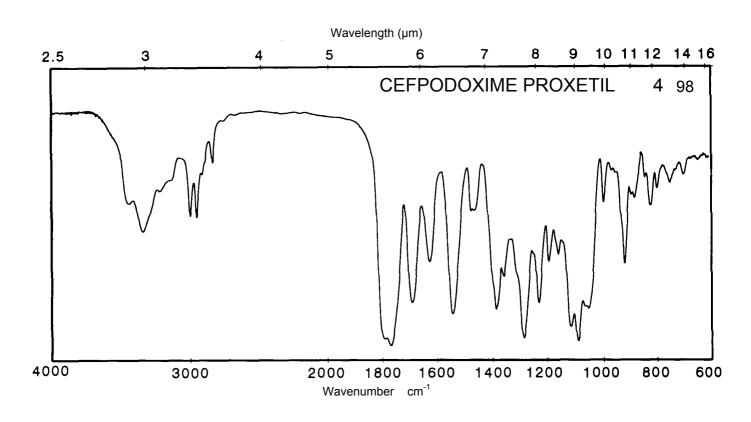





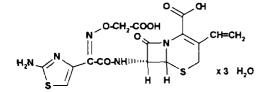


507.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	292 nm 229 nm	Decom- position observed	285 nm	Decom- position observed
E ^{1%} 1cm	438 361		525	
3	22200 18300		26600	

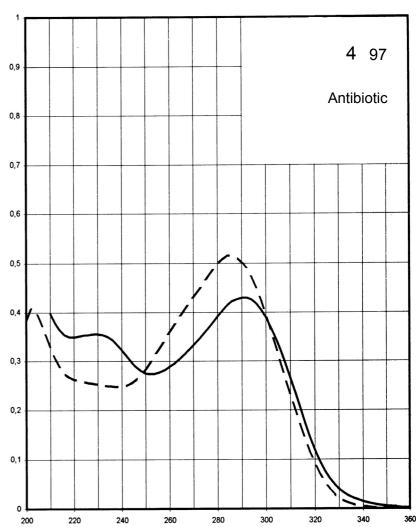

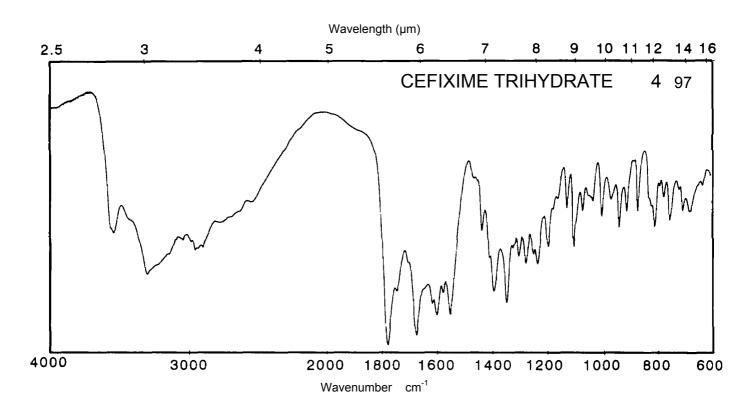





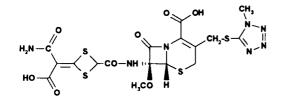


575.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	297 nm	293 nm	298 nm	
E ^{1%} 1cm	429	380	438	
ε	24700	21900	25200	

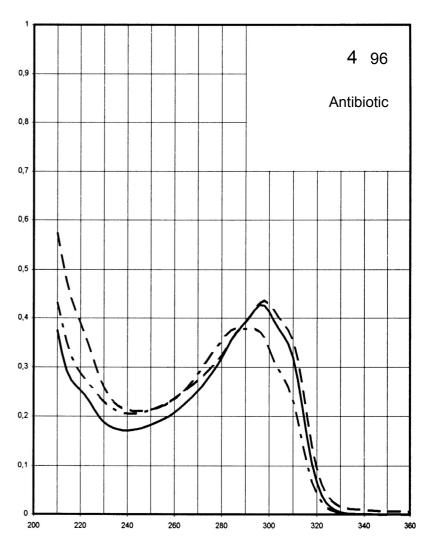

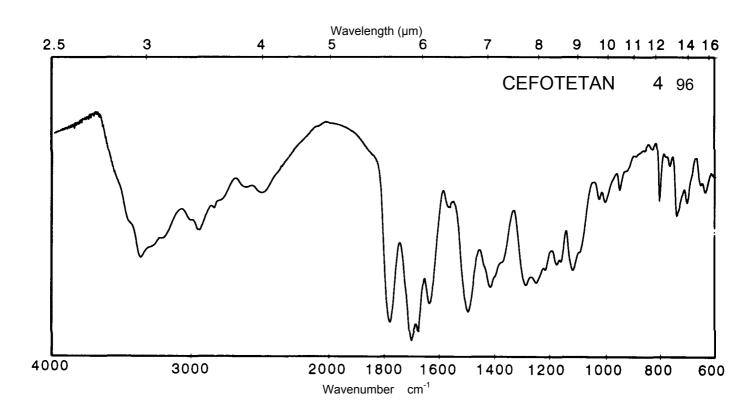




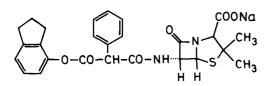



M_r 516.6


Concentration 1 mg / 100 ml 26 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 269 nm	274 nm 268 nm		
E ^{1%} 1cm	32 34	31 34		
3	1640 1740	1610 1760		

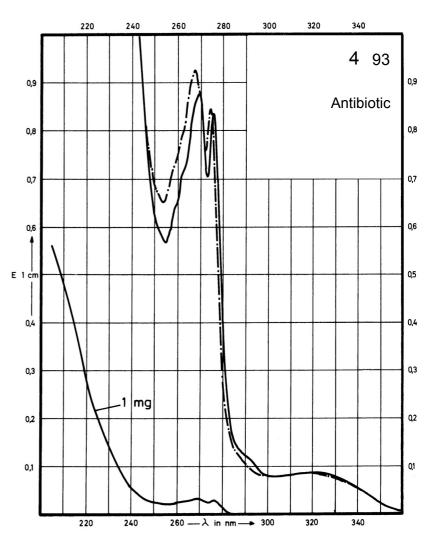

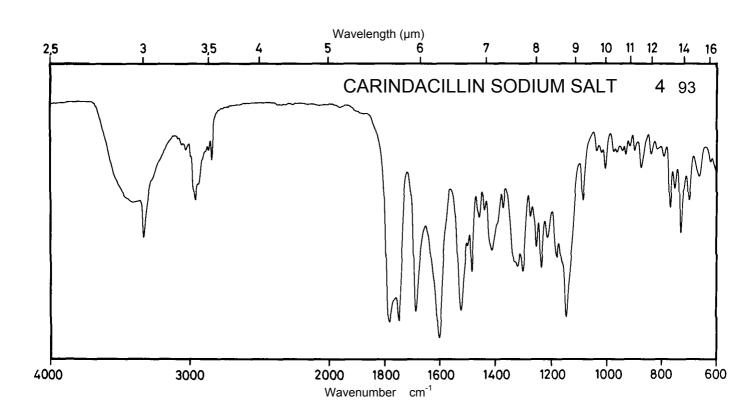

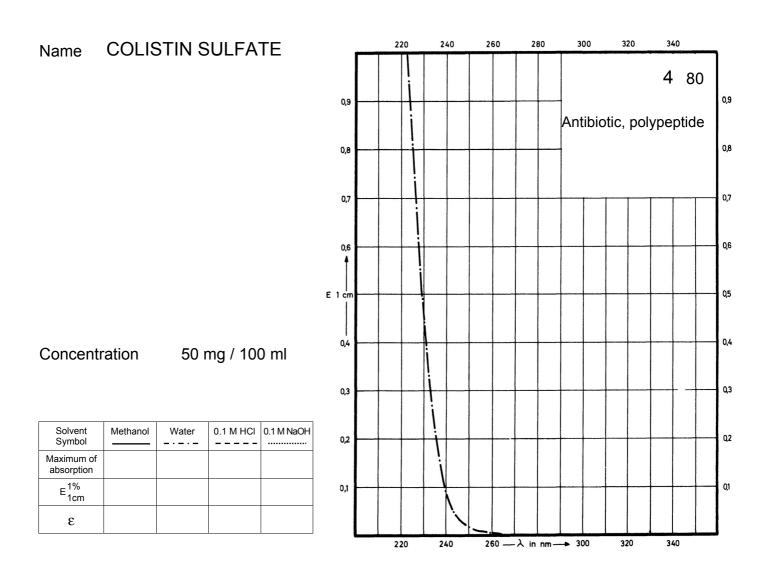

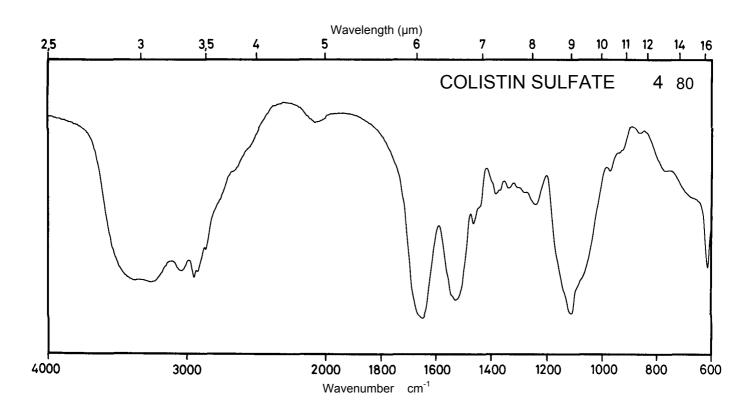

Name CHLORAMPHENICOL SUCCINATE SODIUM SALT

$$\begin{array}{c} \text{COONa} \\ \text{CH}_2 \\ \text{CH-CH-CH}_2 - \text{O-CO-CH}_2 \\ \text{NH-CO-CHCl}_2 \end{array}$$


 M_{r}


445.2

Concentration 2

2.2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		275 nm	280 nm
E ^{1%} _{1cm}	232		215	201
ε	10300		9600	8960

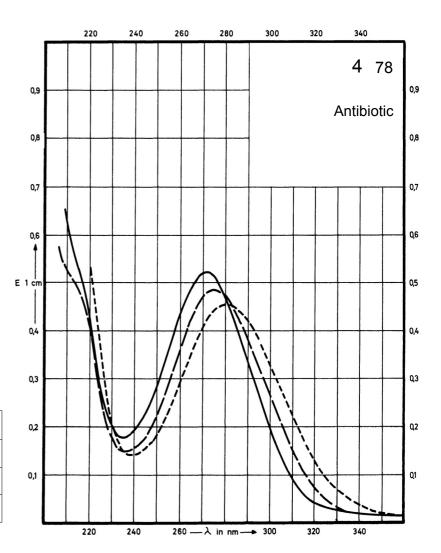

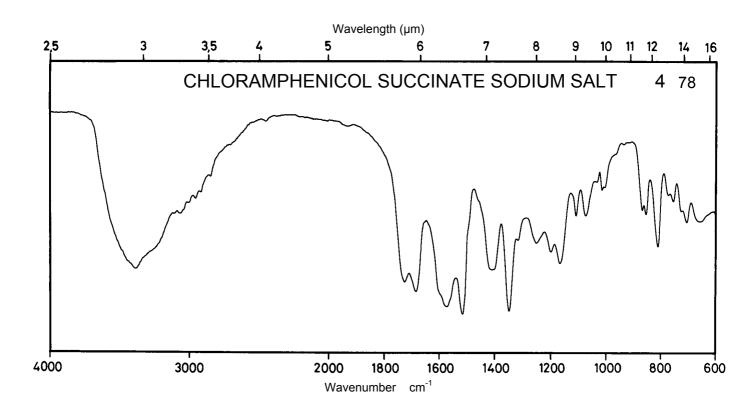

Name CEFTAZIDIME


 M_{r}


636.7

Concentration


2.1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	252 nm	255 nm	259 nm	
E ^{1%} 1cm	395	412	467	
ε	25200	26200	29700	

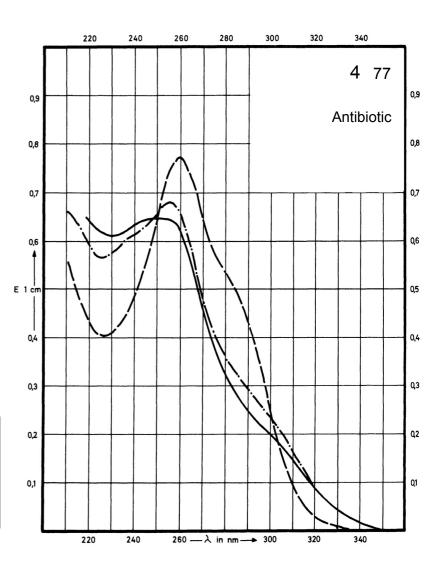

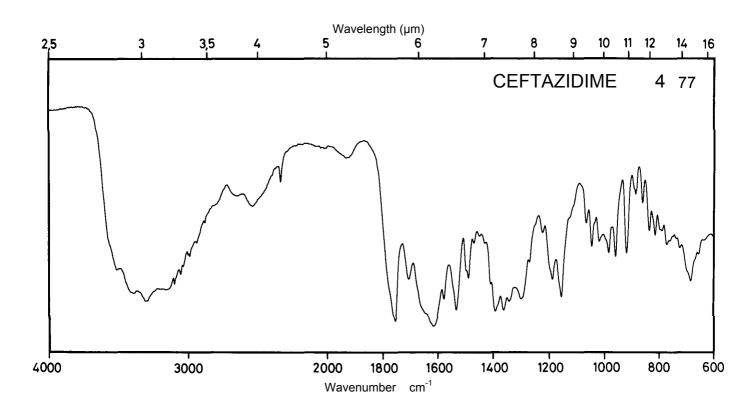





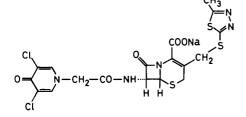


 M_r 570.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm		278 nm	278 nm
E ^{1%} 1cm	521		423	459
3	29730		24150	26200

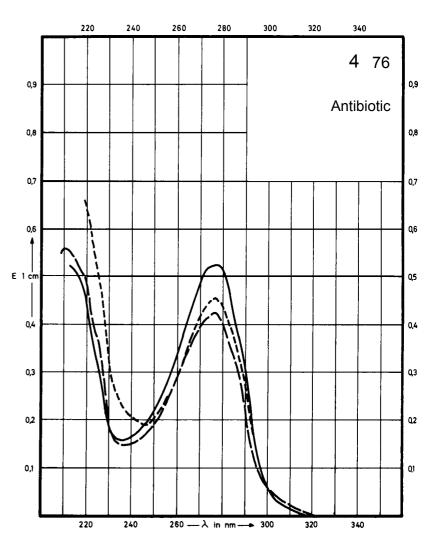

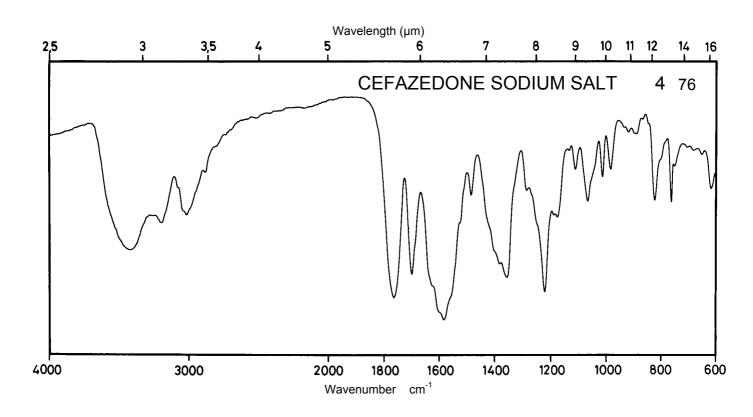


Name CAPREOMYCIN DISULFATE

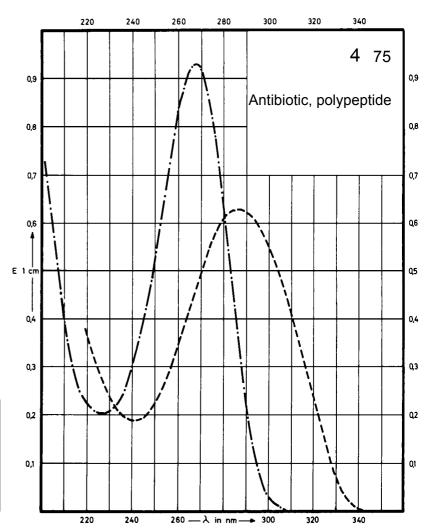


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		268 nm	268 nm	287 nm
E ^{1%} 1cm		306	305	207
ε				

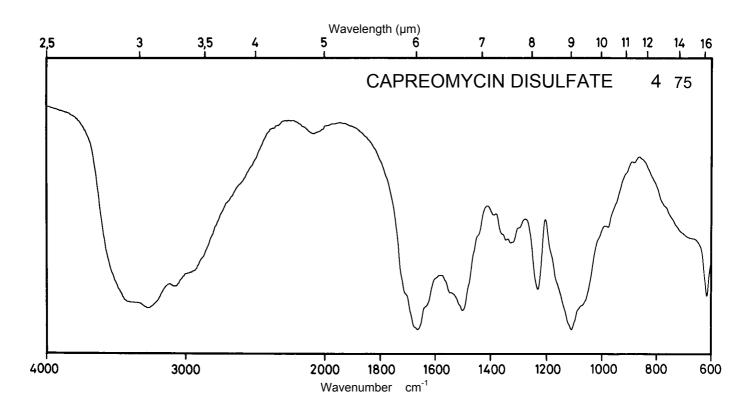





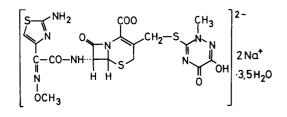

661.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 240 nm		260 nm	273 nm 241 nm
E ^{1%} 1cm	439 529		515	460 478
3	29100 35000		34000	30400 31600

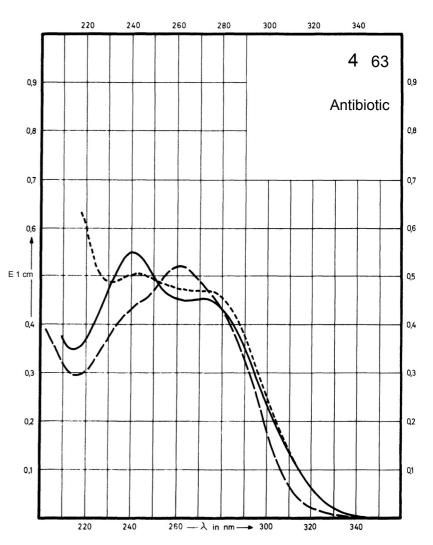


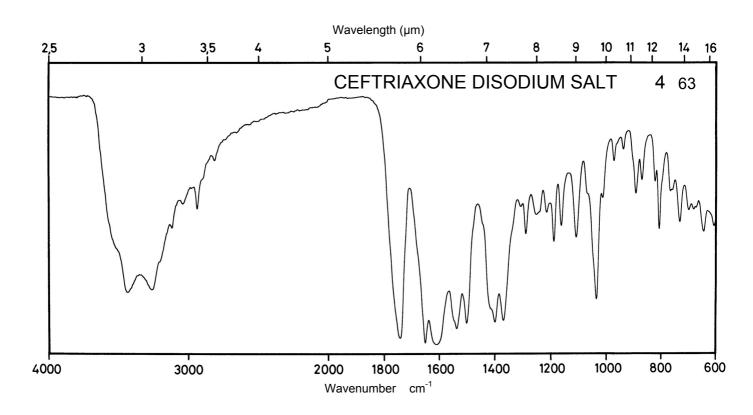





 M_r 401.4

Concentration 4 mg / 100 ml

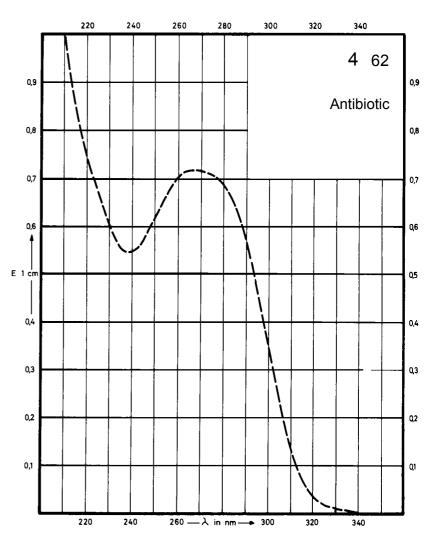

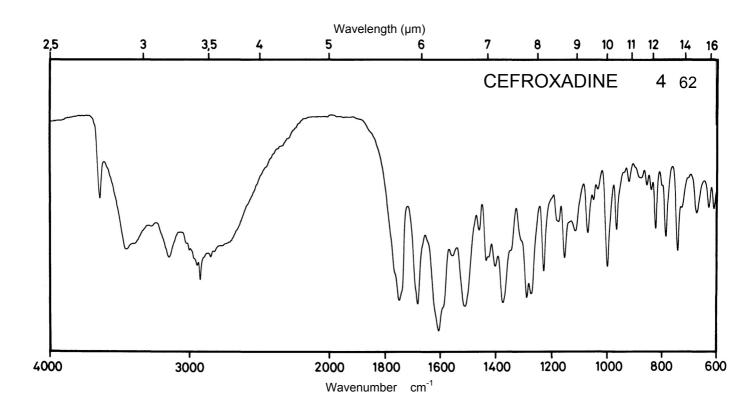





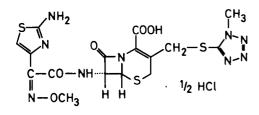

529.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	234 nm		268 nm	Decom- position observed
E ^{1%} 1cm	350		382	
3	18500		20200	

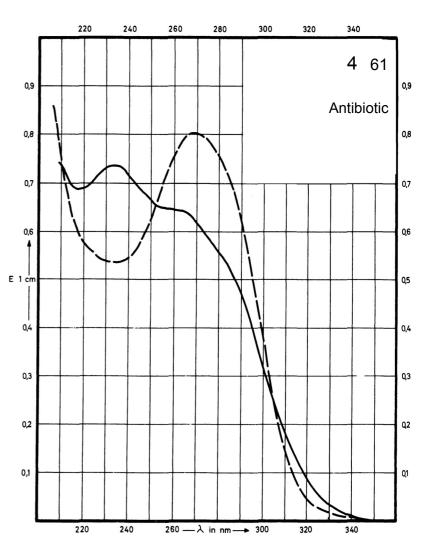


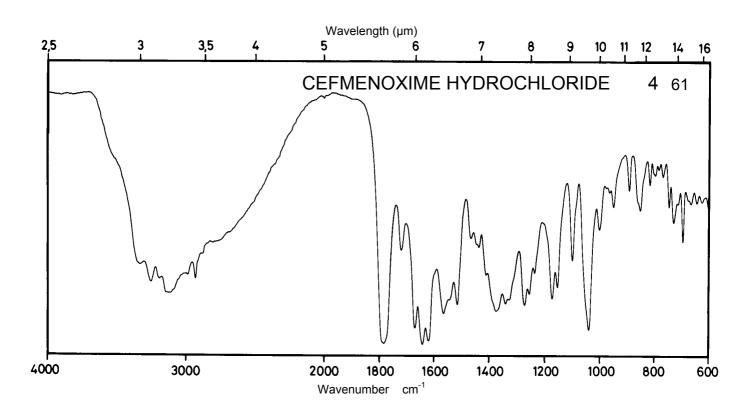





381.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 231 nm		229 nm	249 nm
E ^{1%} 1cm	206 390		350	445
3	7900 14900		13300	17000

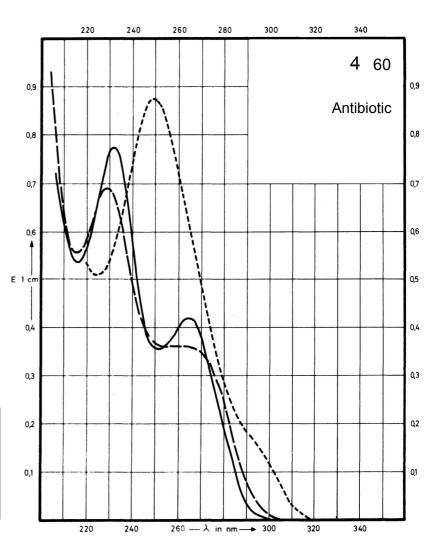

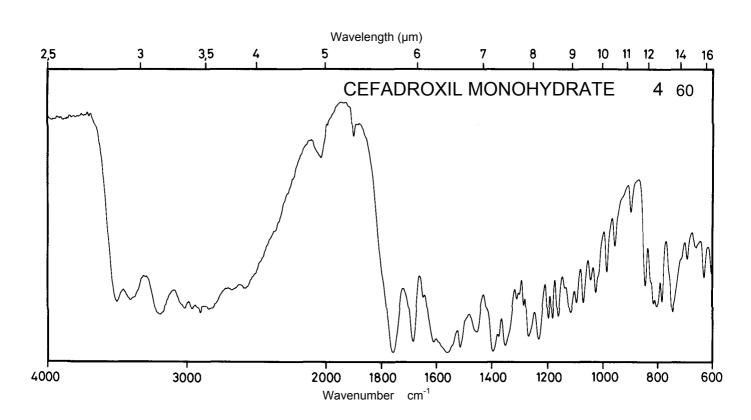





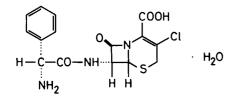


385.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		263 nm	267 nm	
E ^{1%} 1cm		219	198	
3		8470	7620	

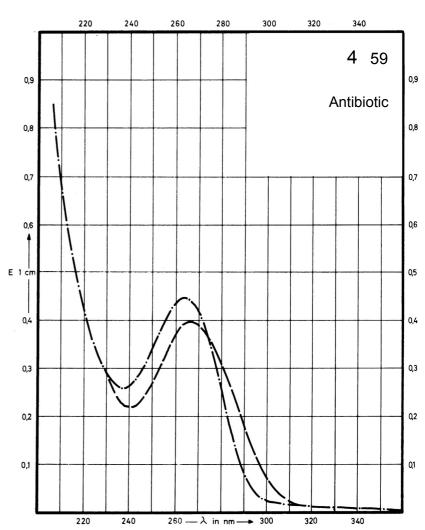

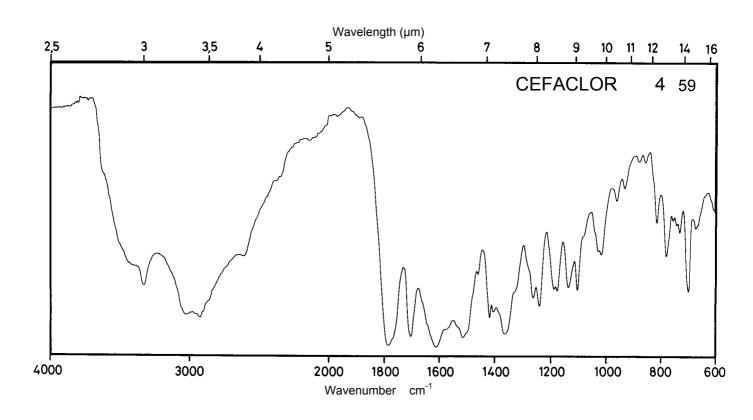





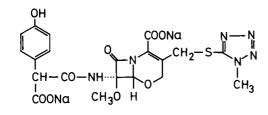


564.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 230 nm	270 nm 227 nm	275 nm 226 nm	Decom- position observed
E ^{1%} 1cm	177 272	169 244	150 207	
ε	10000 15400	9500 13800	8500 11700	

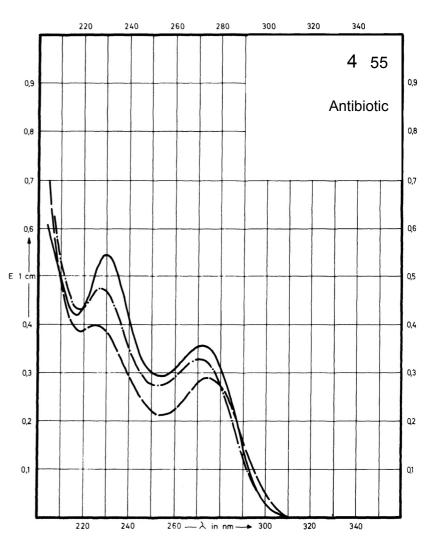

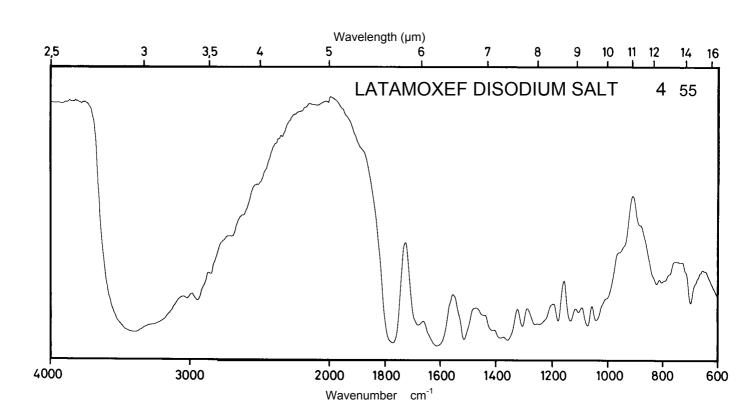





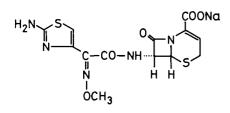


 M_r 405.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	235 nm		260 nm	
E ^{1%} 1cm	Decom- position observed		395	
ε			16000	

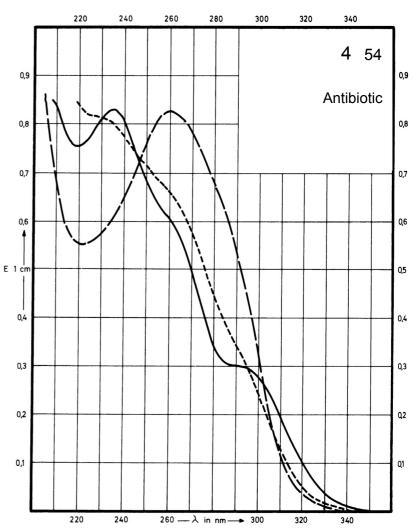

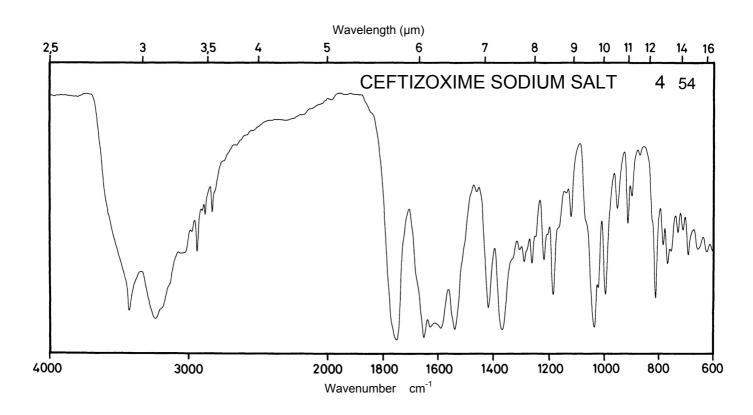





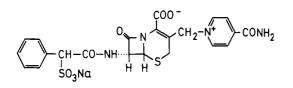



554.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm	263 nm	263 nm	
E ^{1%} 1cm	Decom- position observed	266	261	Decom- position observed
ε		14800	14500	

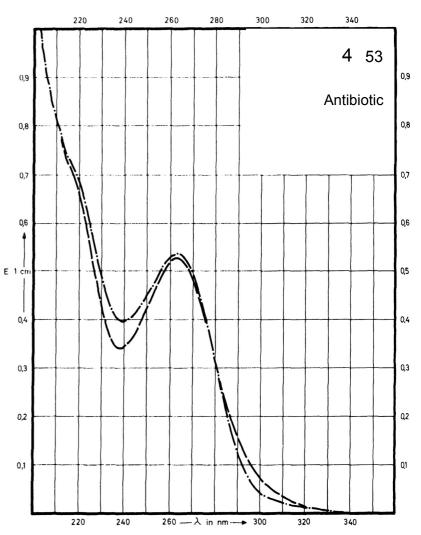


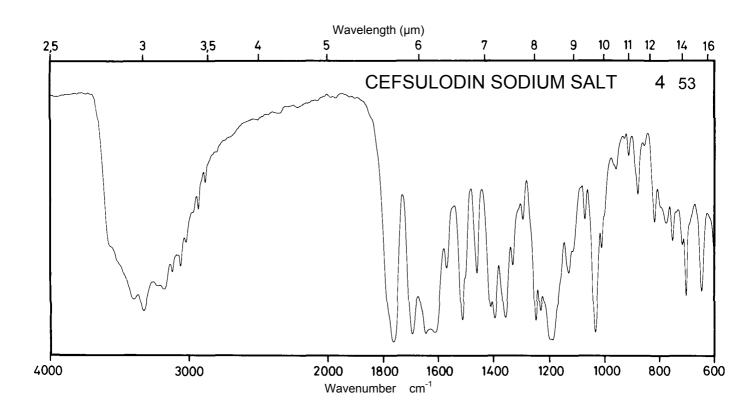

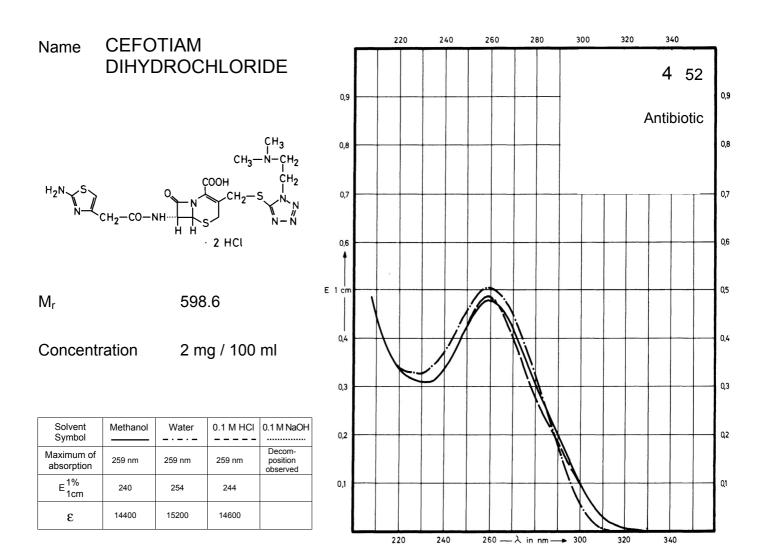

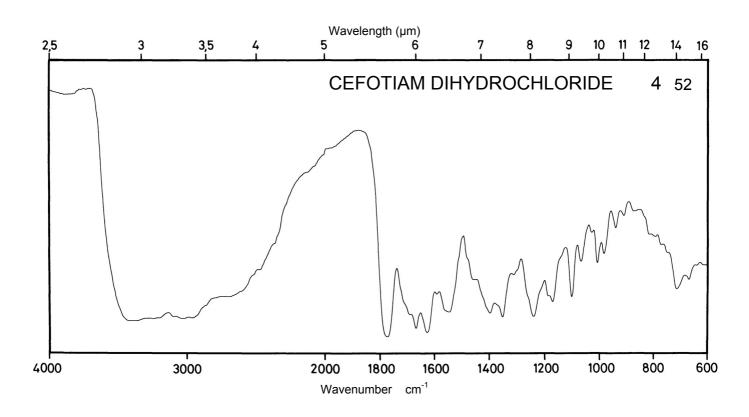





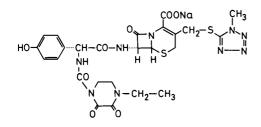


667.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	228 nm	229 nm 266 nm	226 nm	
E ^{1%} 1cm	Decom- position observed	314 218	293	Decom- position observed
3		21000 14600	19600	

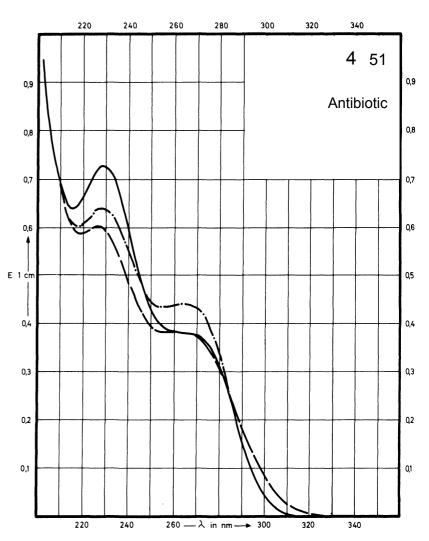


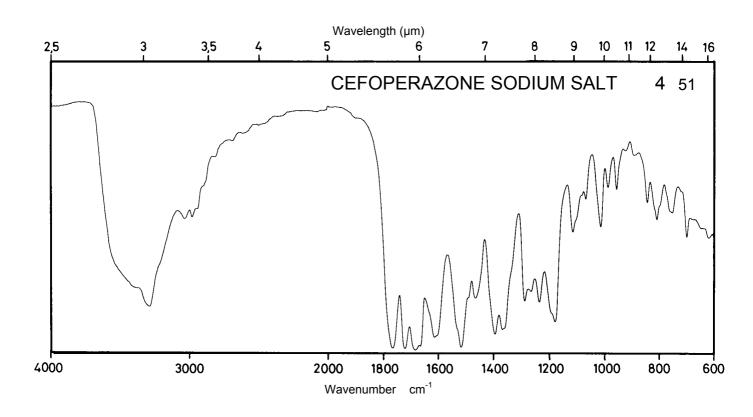




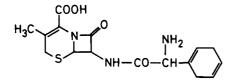

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm		256 nm	Decom- position observed
E ^{1%} 1cm	ca. 140		ca. 124	
ε				

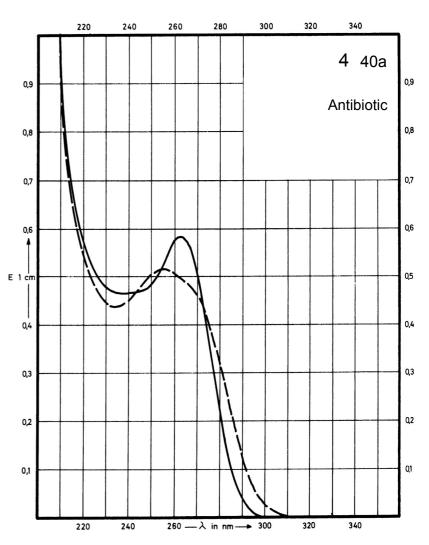

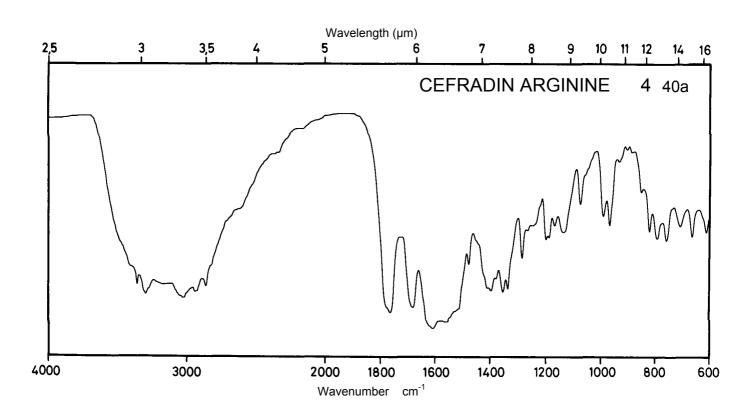



Name CEFRADIN MONO-HYDRATE/ SODIUM CARBONATE

 M_{r}


369.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm		257 nm	Decom- position observed
E ^{1%} 1cm	ca. 190		ca. 170	
3	ca. 7000		ca. 6300	

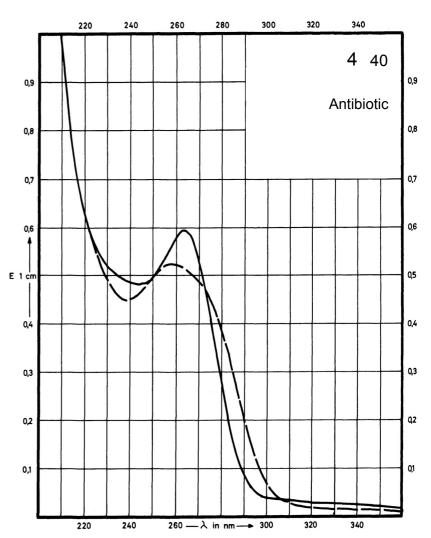

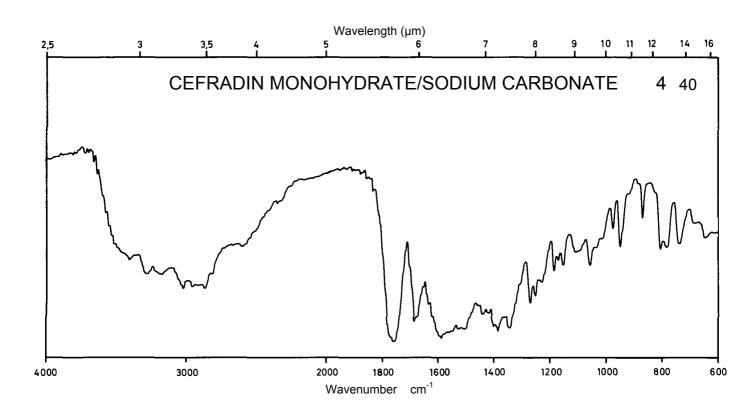





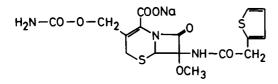


449.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm 237 nm		237 nm	Decom- position observed
E ^{1%} 1cm	180 315		278	
3	8080 14170		12500	

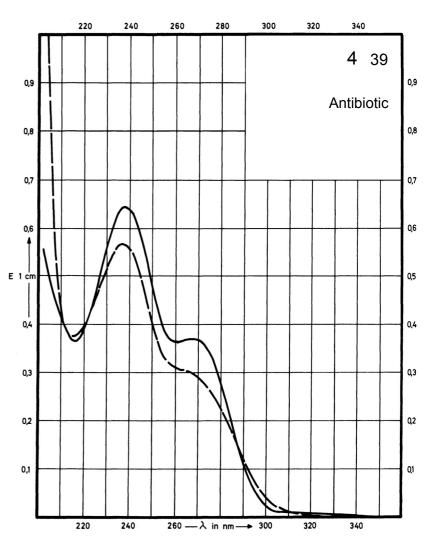

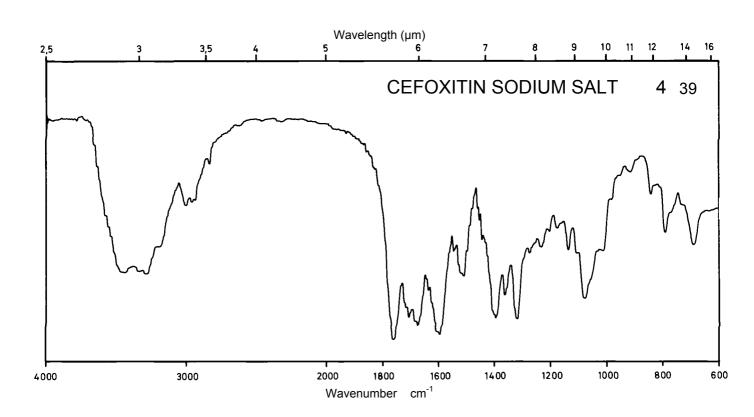





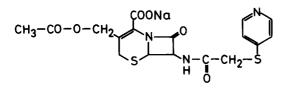


 M_r 445.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm		291 nm 228 nm	Decom- position observed
E ^{1%} 1cm	Decom- position observed		420 282	
3			18700 12570	

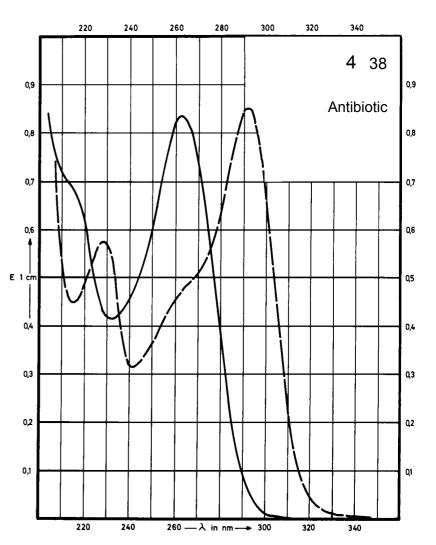

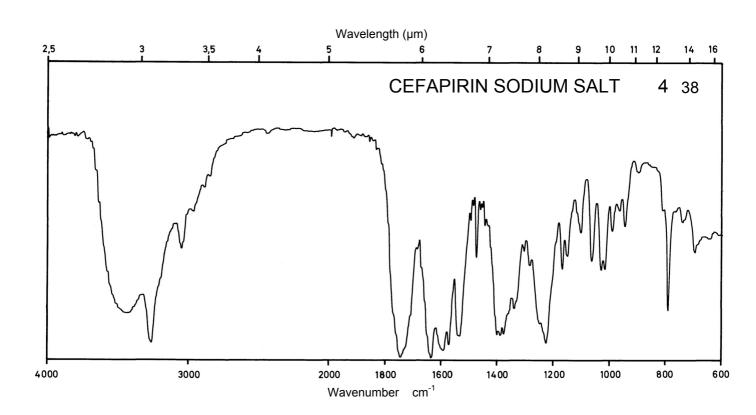




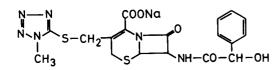


 M_r


484.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm		270 nm	Decom- position observed
E ^{1%} 1cm	Decom- position observed		150	
ε			7250	

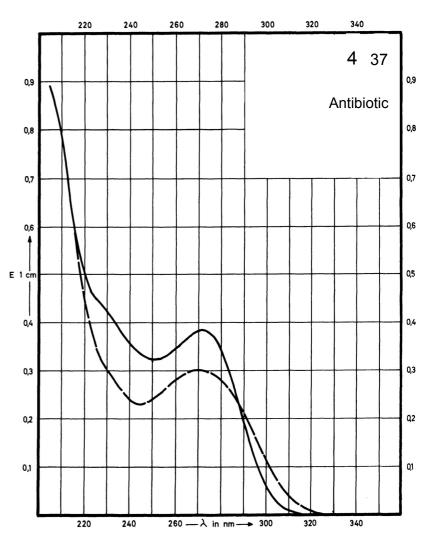

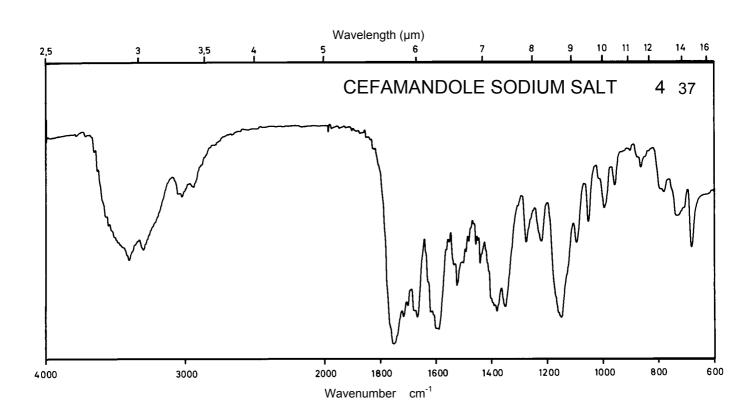





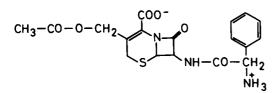


 M_r 405.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		262 nm	258 nm	Decom- position observed
E ^{1%} 1cm		210	194	
3		8510	7860	

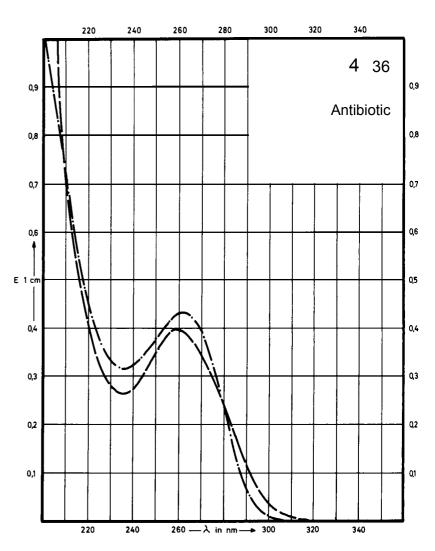

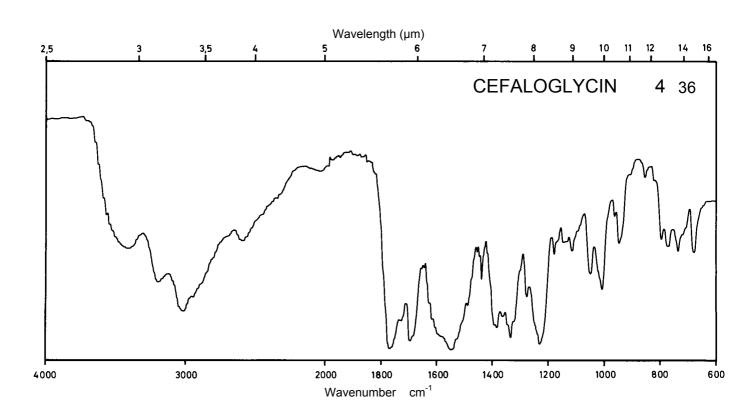



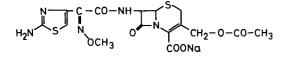



477.5

 M_r


Concentration 2 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 235 nm 235 nm 262 nm absorption E^{1%} 1cm 397 370 391 17640 ε 18940 18680

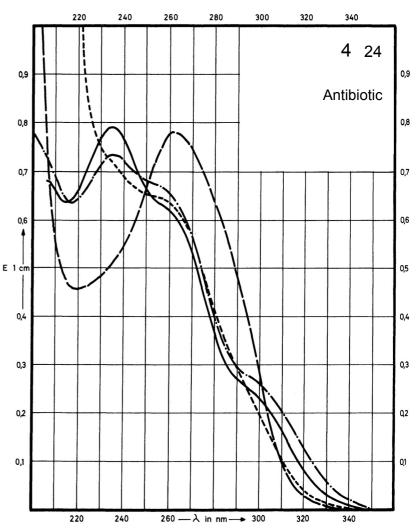

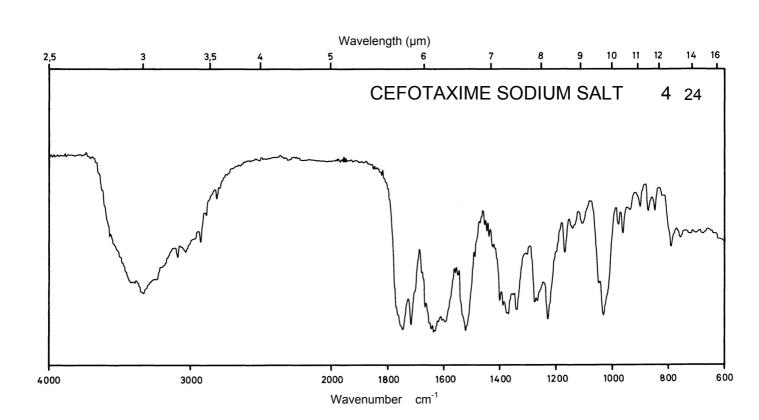




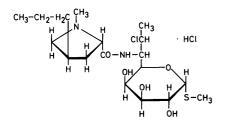


 M_r


531.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

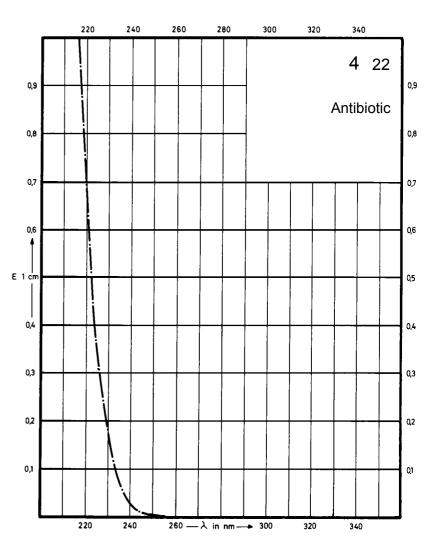

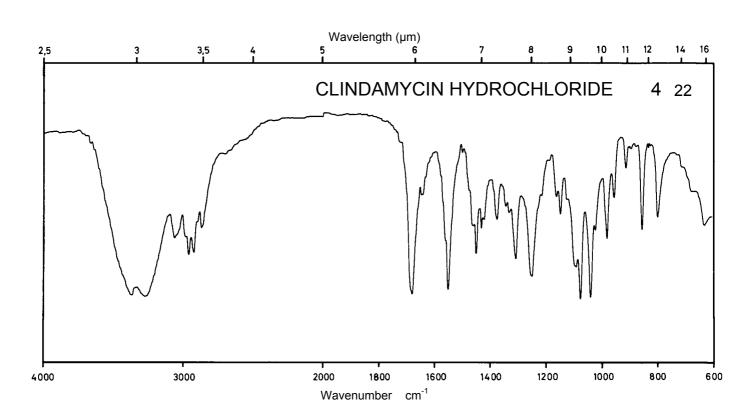





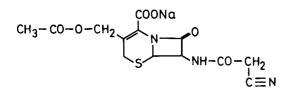


361.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	260 nm	260 nm	Decom- position observed
E ^{1%} 1cm		243	222	
ε		8770	8000	

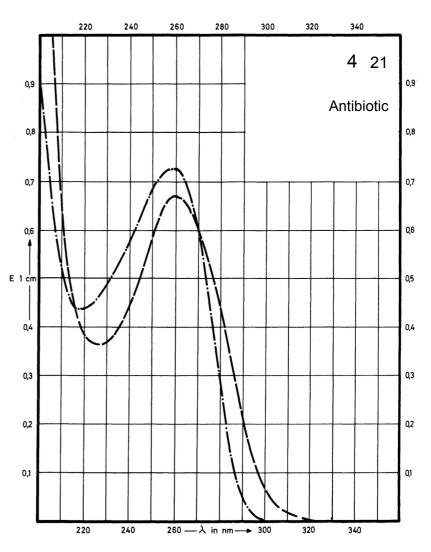

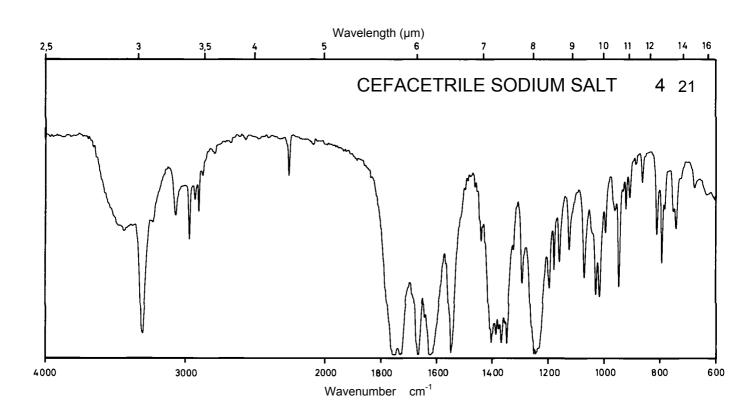





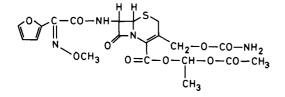


510.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm	280 nm	280 nm	Decom- position observed
E ^{1%} _{1cm}	400	376	376	
ε	20400	19200	19200	

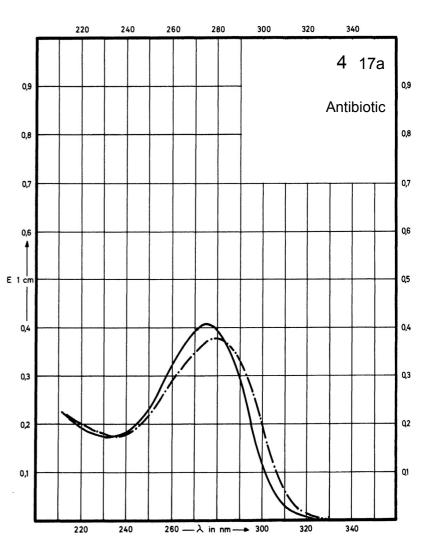

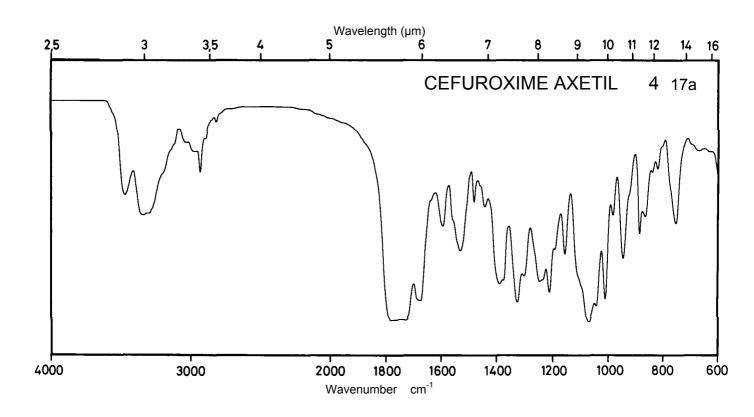





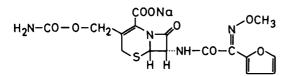


 M_r 446.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm		280 nm	272 nm
E ^{1%} _{1cm}	410		400	440
ε	18300		17860	19640

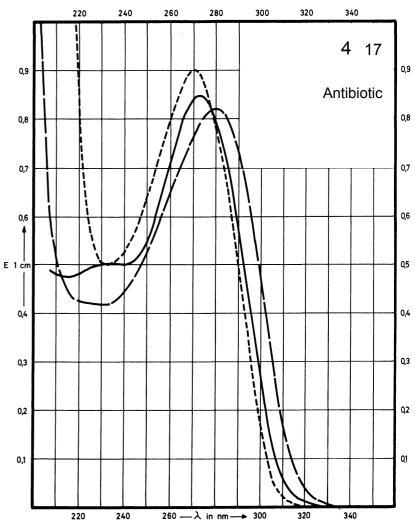

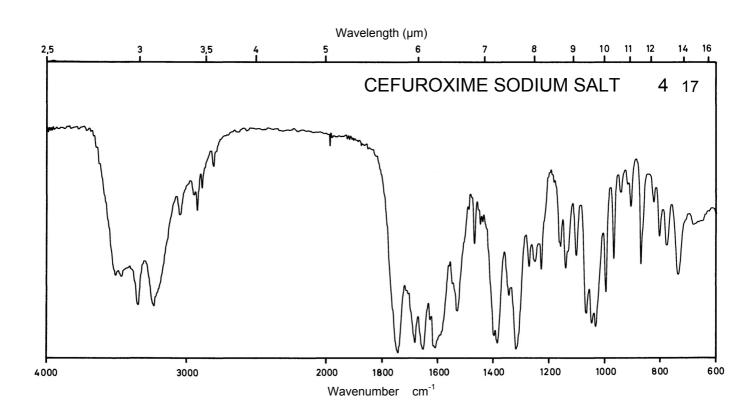





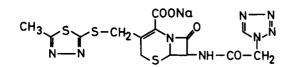


 M_r 476.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			266 nm	
E ^{1%} 1cm			214	
ε			10200	

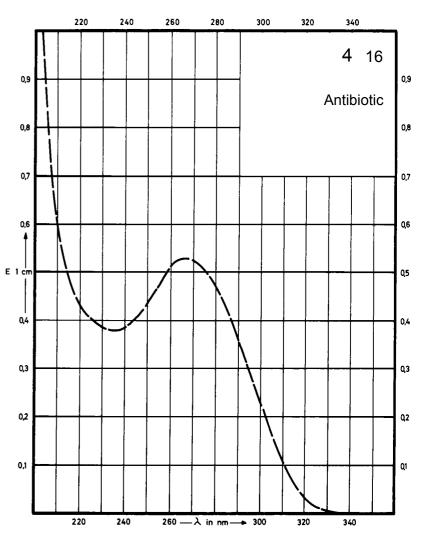

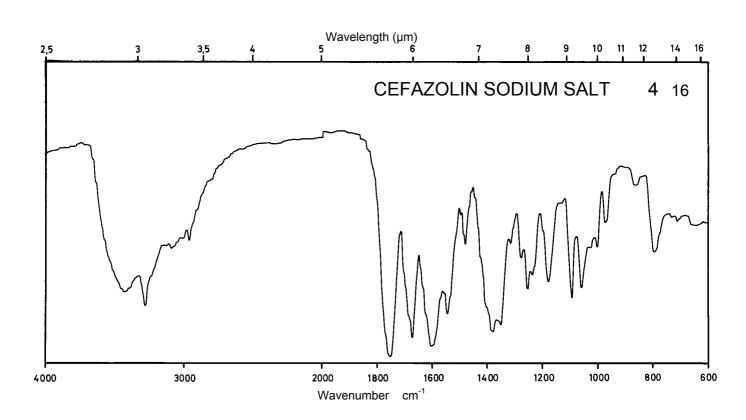





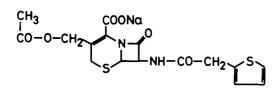

 M_r 418.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 236 nm	236 nm		
E ^{1%} 1cm	220 371	339		
ε	9200 15520	14180		

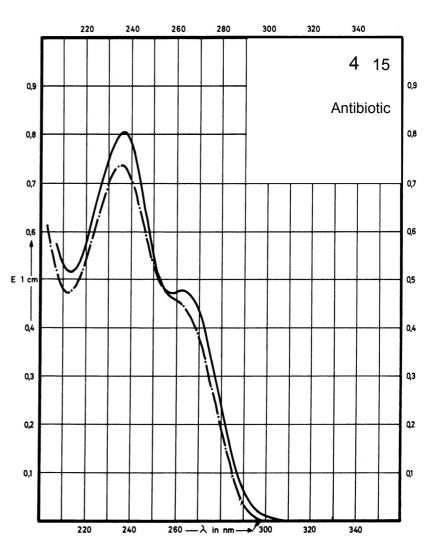

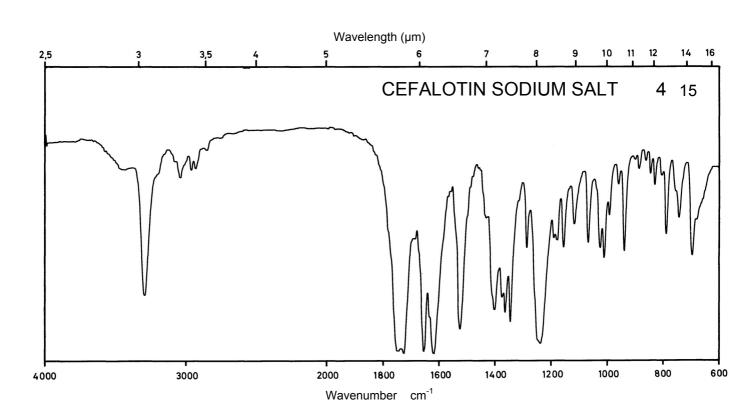



Name CHLORTETRACYCLINE HYDROCHLORIDE

 M_{r}


515.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm 232 nm		266 nm 228 nm	284 nm 254 nm
E ^{1%} 1cm	349 348		360 338	280 287
3	17990 17940		18550 17420	14430 14790

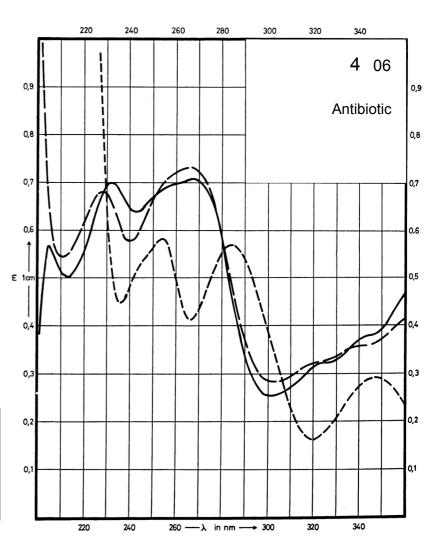

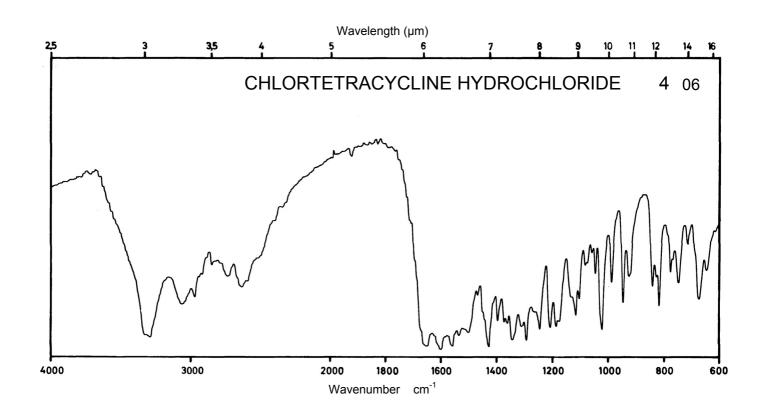



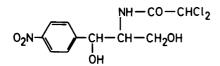



323.1

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm		277 nm	277 nm
E ^{1%} 1cm	297		289	291
ε	9600		9340	9400

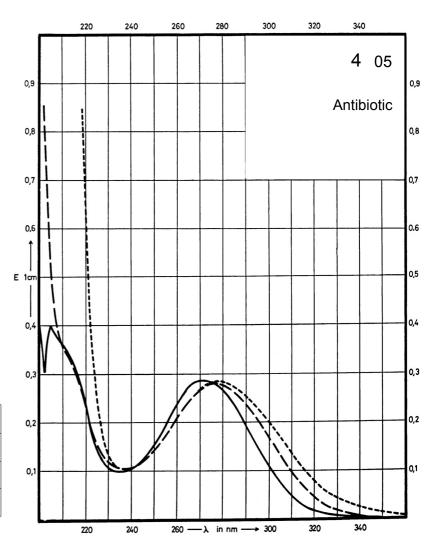

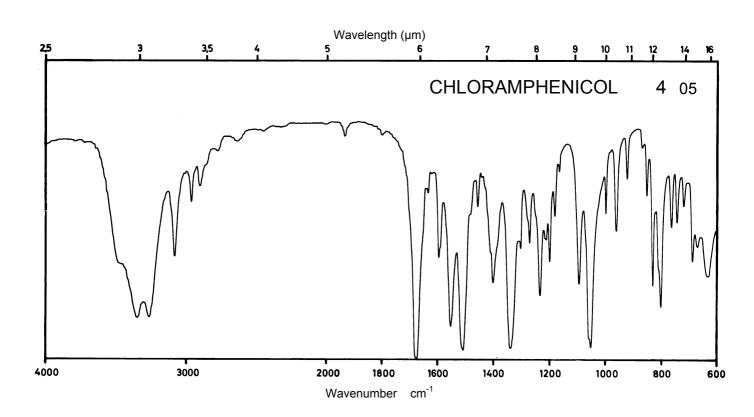



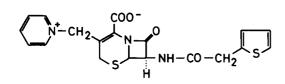



 M_r

415.5


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	233 nm		254 nm 240 nm	261 nm
E ^{1%} 1cm	374		327 335	88
ε	15540		13590 13920	3660

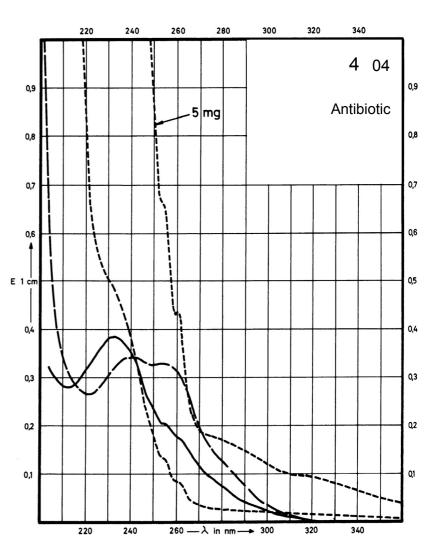

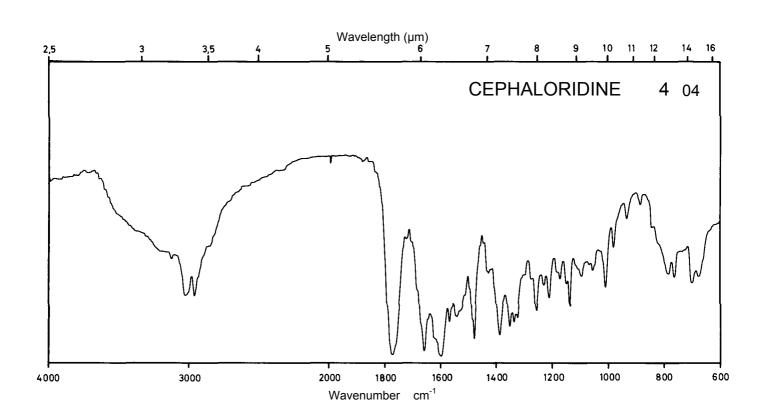




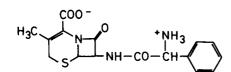


347.4


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm	260 nm		
E ^{1%} 1cm	226	219		
ε	7850	7610		

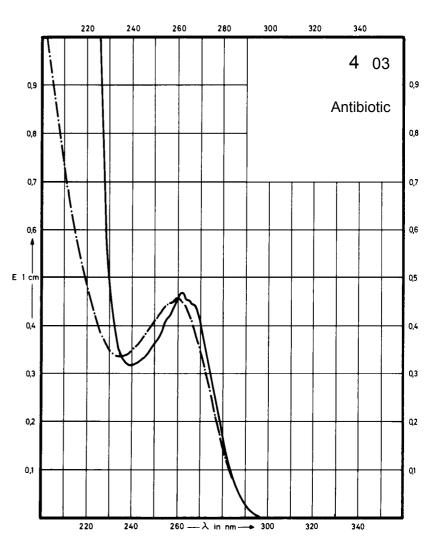

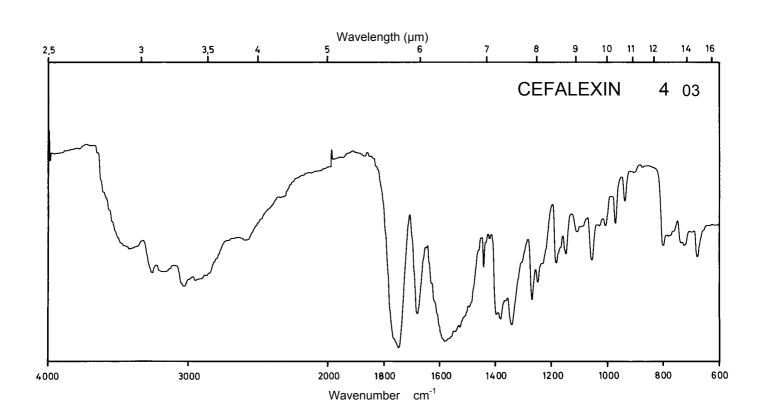




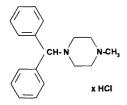


302.9


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm 254 nm 225 nm	225 nm	225 nm	
E ^{1%} 1cm	15 14 341	359	393	
ε	462 435 10300	10900	11900	

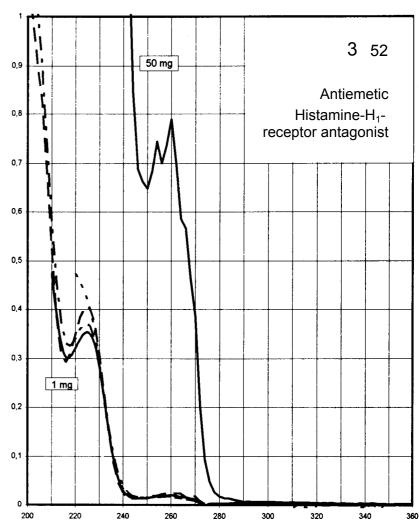

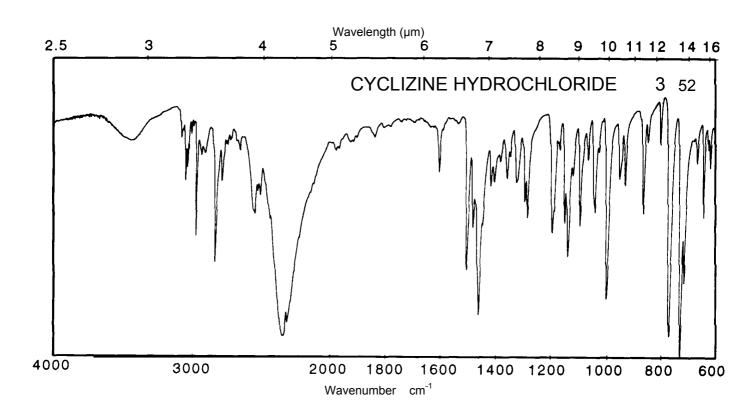



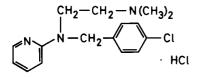


 M_r


Concentration

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 306 nm 244 nm 314 nm 238 nm 312 nm 247 nm absorption E^{1%} 1cm 148 497 257 478 142 485 4820 16200 8400 15600 4620 15800 ε


326.3

1 mg / 100 ml

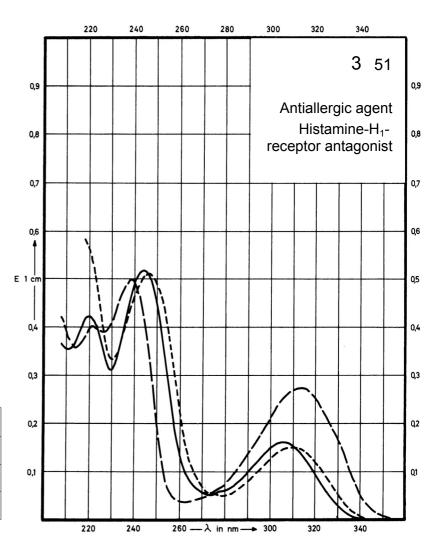

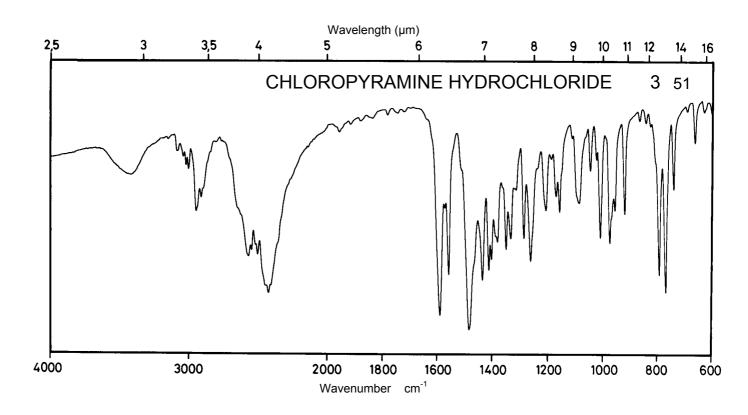

Name CLEMIZOLE HYDROCHLORIDE

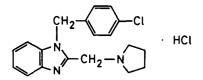



 M_{r}

362.3


Concentration 2.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm 275 nm 251 nm		275 nm	282 nm 274 nm 253 nm
E ^{1%} 1cm	158 176 220		279	160 182 208
3	5700 6360 7960		10100	5800 6580 7540

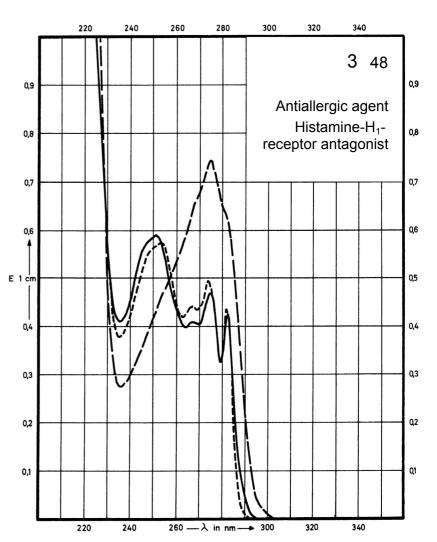

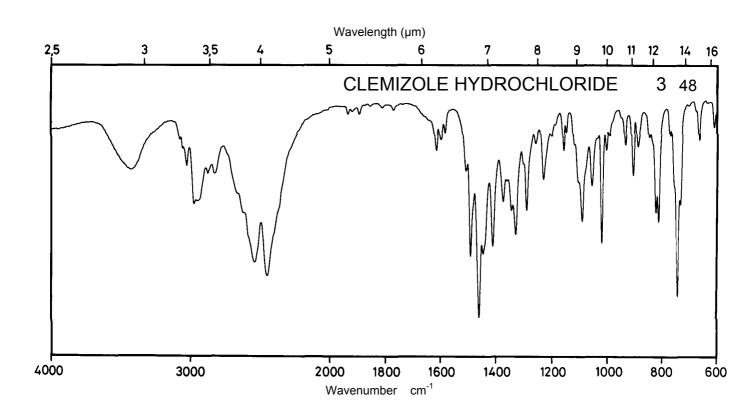



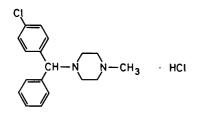



 M_r

337.3

Concentration 1.6 mg / 100 ml 2 mg / 100 ml 40 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm 229 nm		231 nm	229 nm
E ^{1%} 1cm	16.8 442		505	458
ε	570 14900		17000	15500

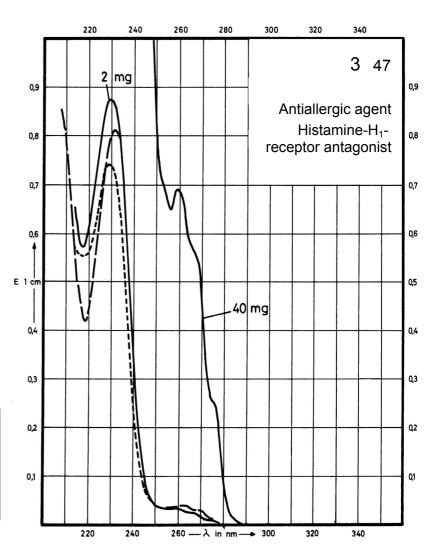

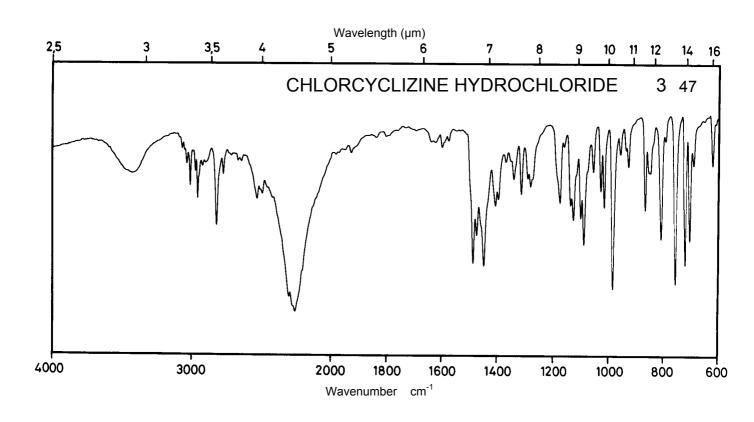


Name 8-CHLORTHEO-PHYLLINE


 M_{r}

214.6


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm		276 nm	278 nm
E ^{1%} 1cm	541		576	629
ε	11600		12400	13500

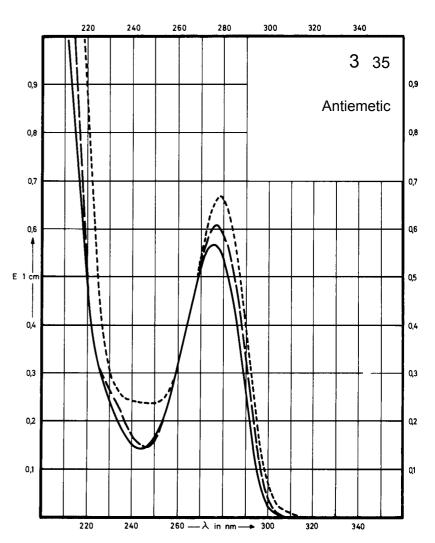

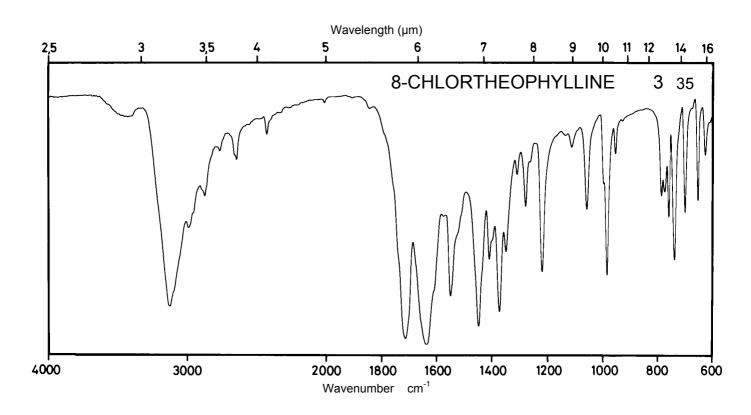




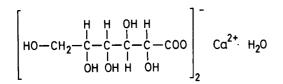


 M_{r} 448.4


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

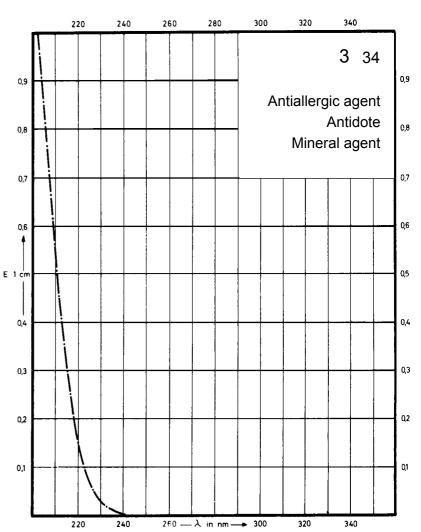

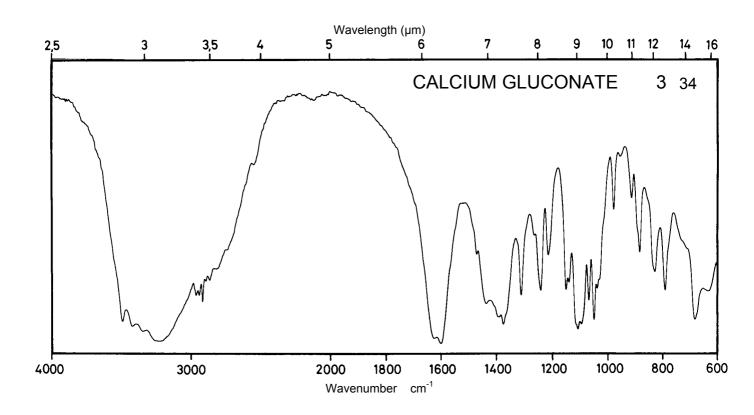

Name CHLORPHENOXAMINE HYDROCHLORIDE

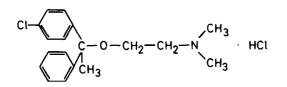



 M_r

340.3


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 259 nm		259 nm	
E ^{1%} 1cm	6.1 13.6		14.7	
ε	209 464		500	

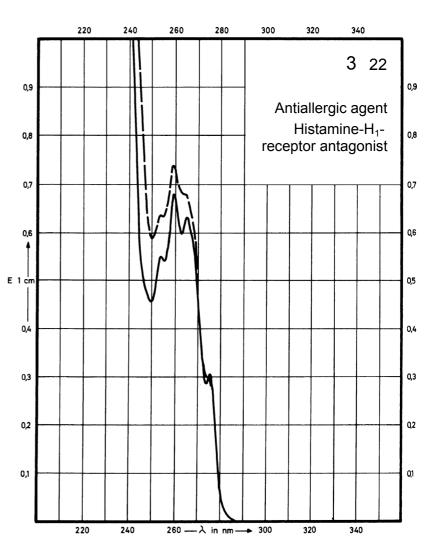

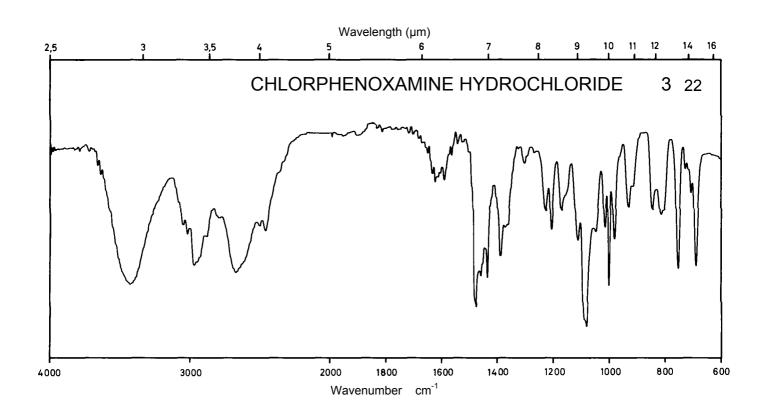

Name CHLORPHENOXAMINE

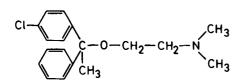



 M_{r}

303.8

Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 259 nm		259 nm	
E ^{1%} 1cm	7.5 16.0		16.9	
3	230 490		510	

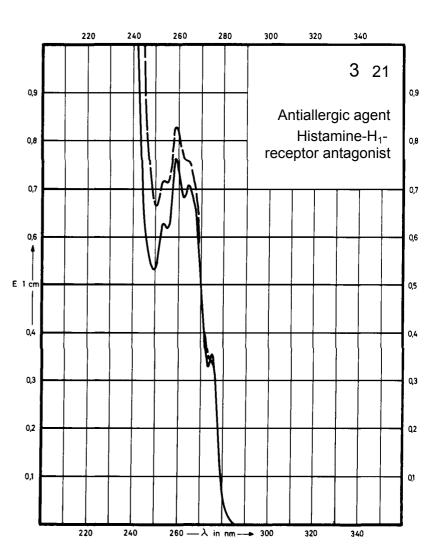

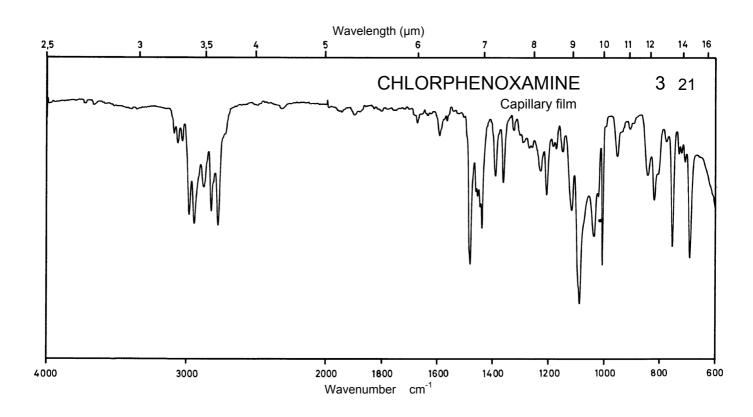





460.0


Concentration 30 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm			
E ^{1%} 1cm	20.3			
ε	935			

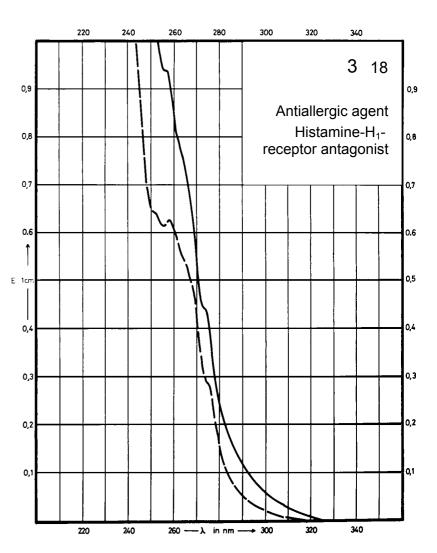

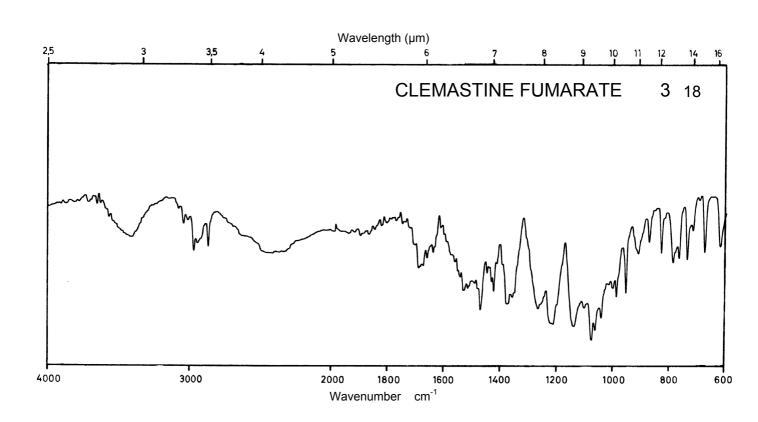


Name CYPROHEPTADINE HYDROCHLORIDE

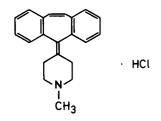


 M_{r} 323.9

Concentration 0.4 mg / 100 ml 0.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm 223 nm		284 nm 223 nm	284 nm
E ^{1%} 1cm	410 1345		430 1390	377
ε	13280 43560		13930 45020	12210

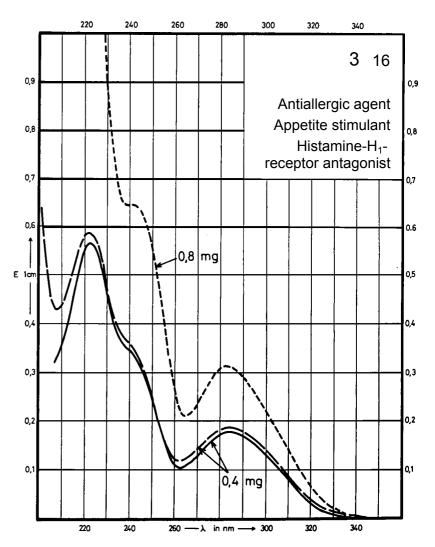

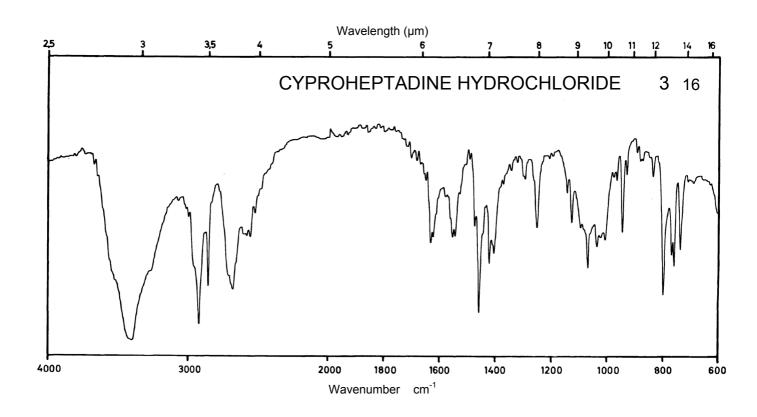


Name CHLORPHENAMINE MALEATE


 M_{r}

390.9


Concentration 2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm		264 nm	262 nm
E ^{1%} 1cm	140		218	144
ε	5470		8520	5630

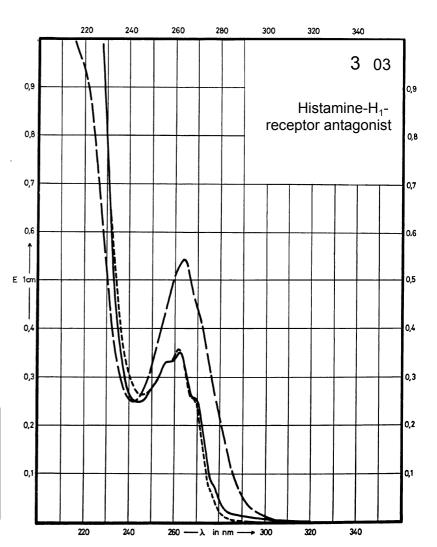

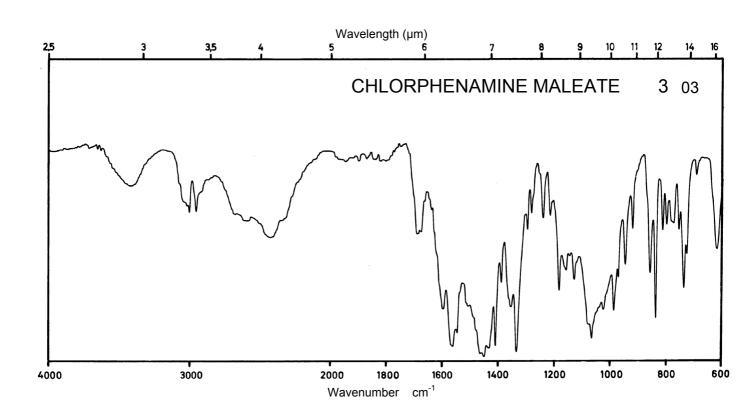




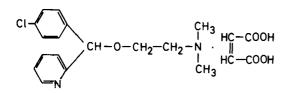


 M_r 406.9

Concentration 2.1 mg / 100ml 5.3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm		263 nm	260 nm
E ^{1%} 1cm	136		247	134
ε	5530		10050	5450

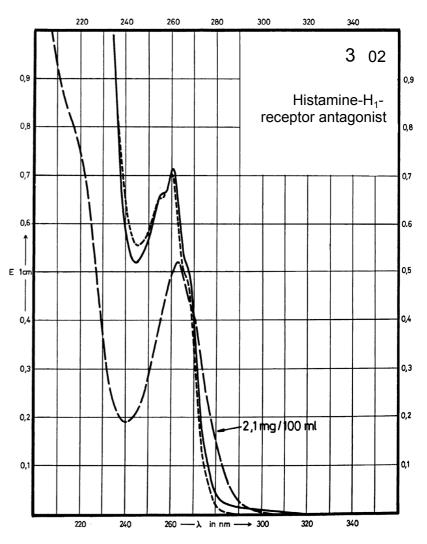

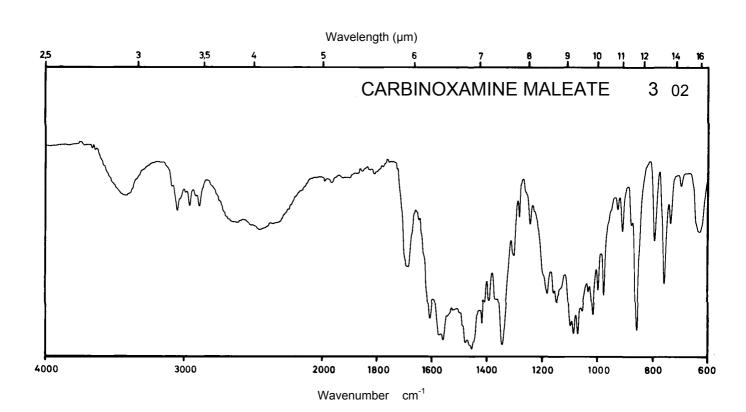


Name CLIOXANIDE


 M_{r}

541.5

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	232 nm			362 nm 282 nm 240 nm
E ^{1%} 1cm	602			172 270 561
3	32600			9310 14620 30380

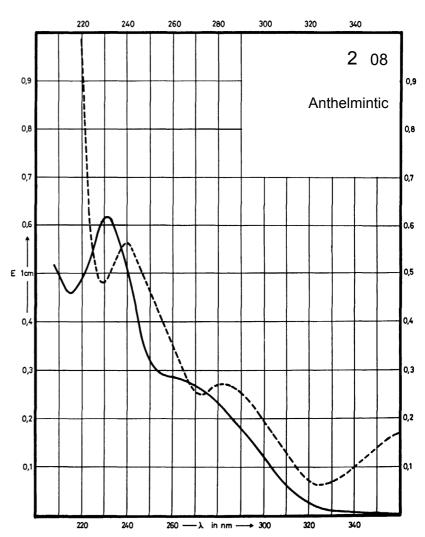

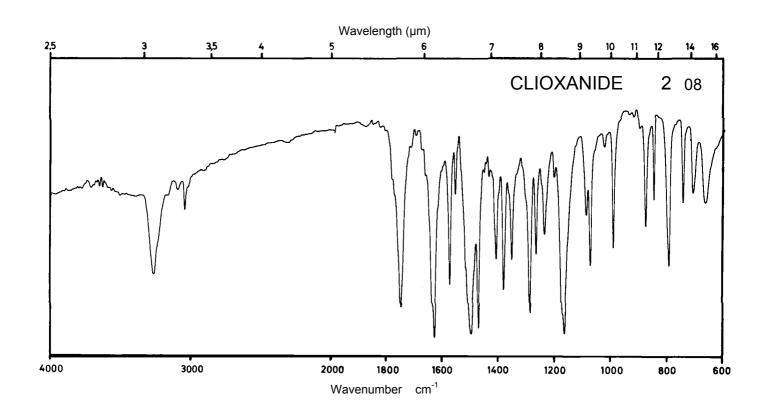


Name CAMBENDAZOLE


 M_{r}

302.4


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	317 nm 238 nm		319 nm 235 nm	314 nm 237 nm
E ^{1%} 1cm	702 813		689 620	649 936
ε	21230 24590		20840 18750	19630 28300

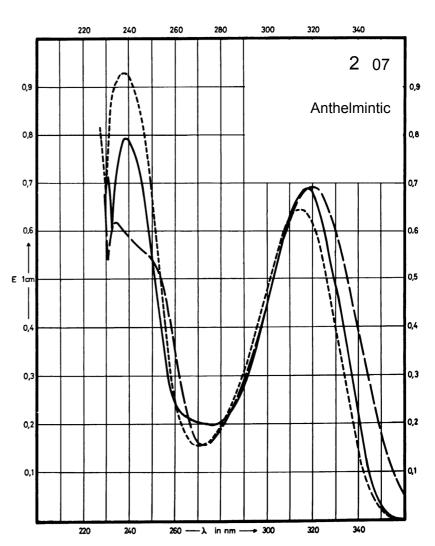

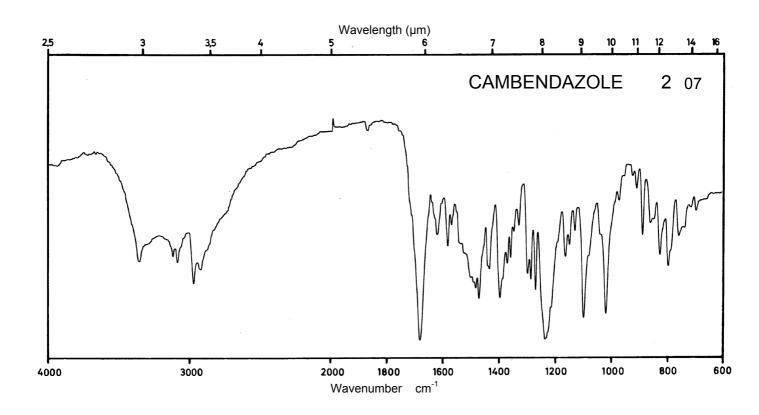




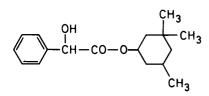


276.4


Concentration 60 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm			
E ^{1%} 1cm	6.05 7.72 6.40			
3	167 213			

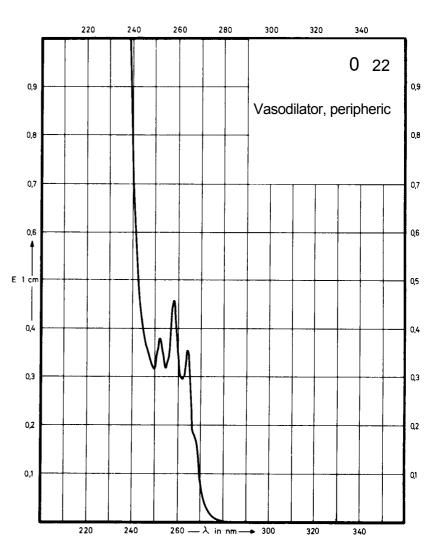

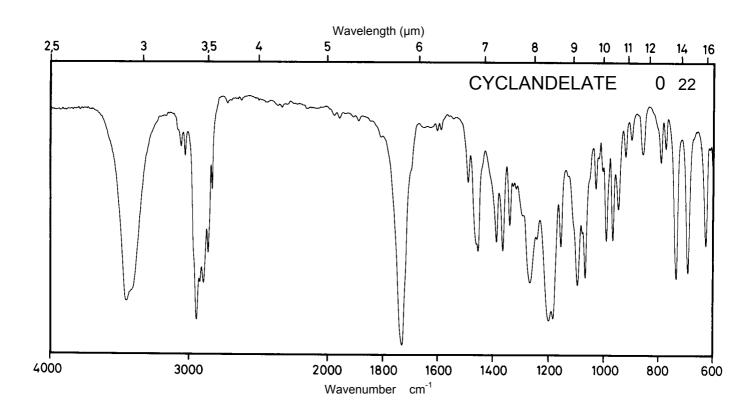




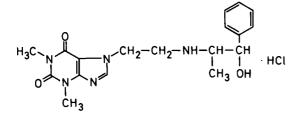


393.9


2 mg / 100 ml Concentration

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm		274 nm	274 nm
E ^{1%} 1cm	222		232	203
ε	8800		9100	8000

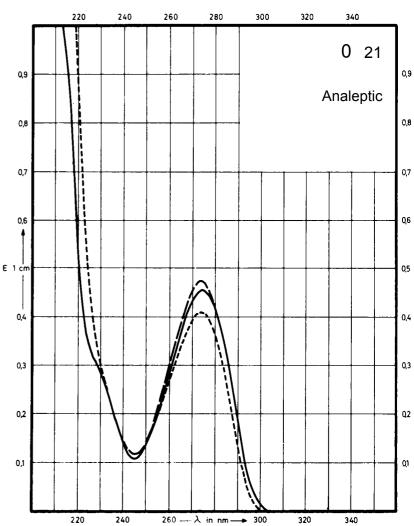

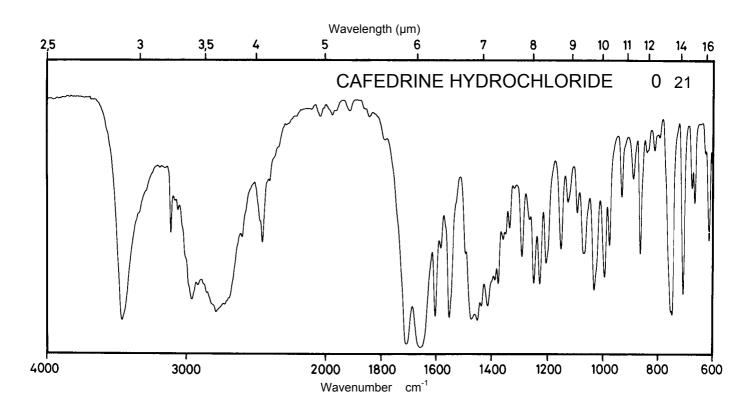





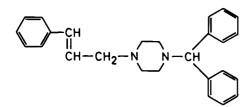


 M_r 368.5

Concentration 1 mg / 100 ml 20 mg / 100 ml

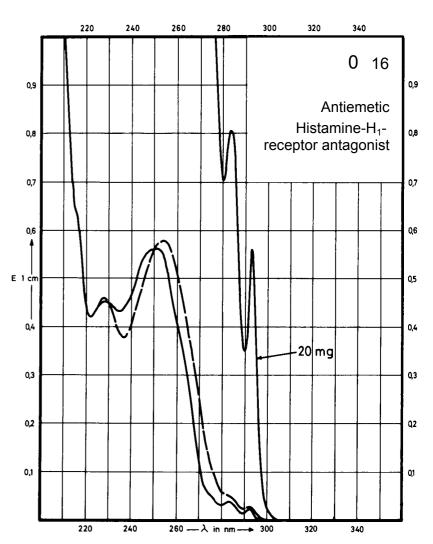

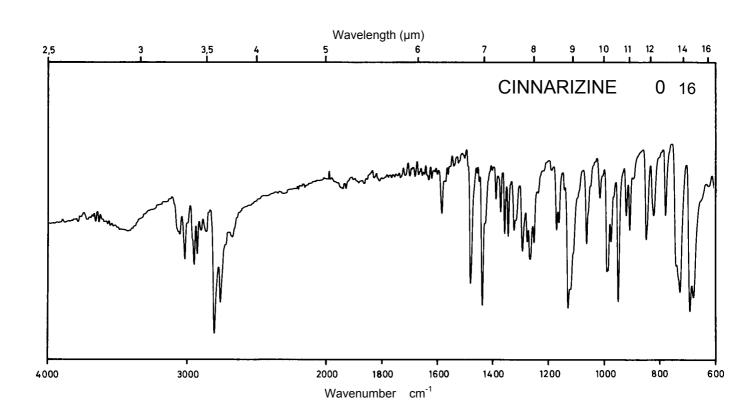

Name CAMPHOR

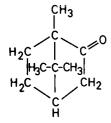


 M_{r}


152.2

Concentration


240 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm			
E ^{1%} 1cm	2.07			
3	31.5			

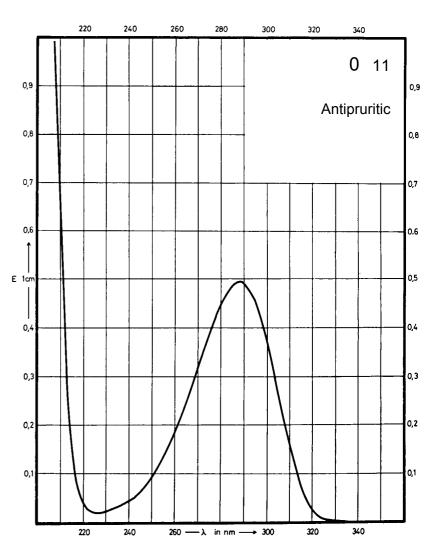

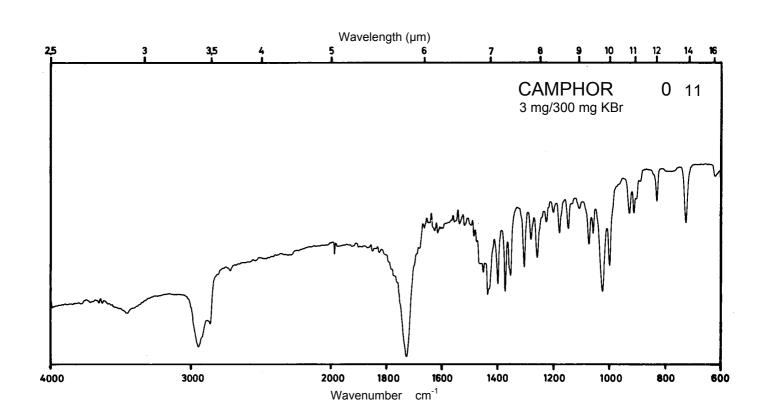




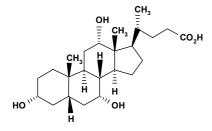


 M_r 408.6

Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

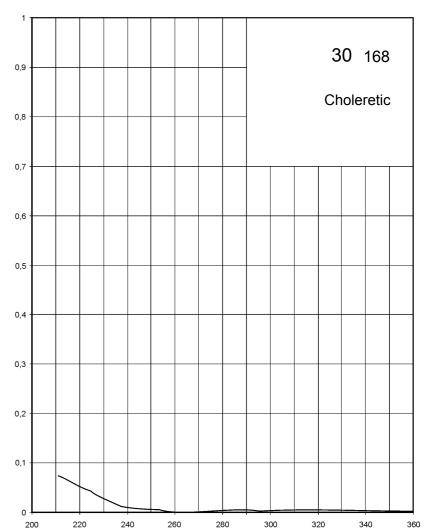

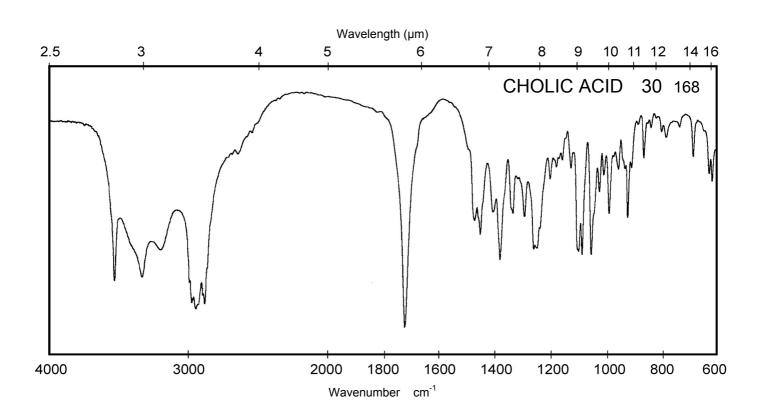




C19H20CIN3 · C16H18N2O4S


 M_r

Concentration 2.5 mg / 100 ml


660.2

Ī	Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
ŀ		283 nm	281 nm		
	Maximum of absorption	275 nm 252 nm	274 nm 252 nm		
ŀ	10/	94	88		
	E ^{1%} 1cm	96	107		
L	ICIII	128	122		
		6200	5800		
	3	6360	7040		
	-	8470	8050		

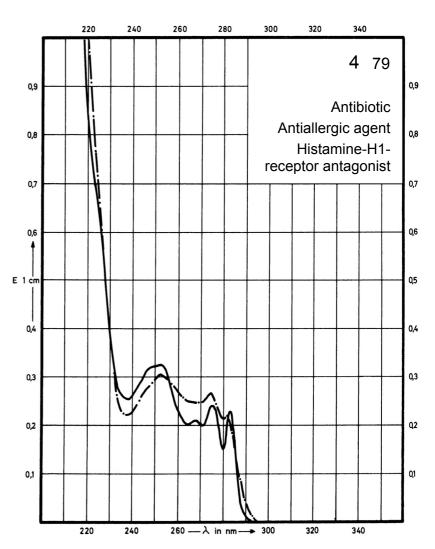

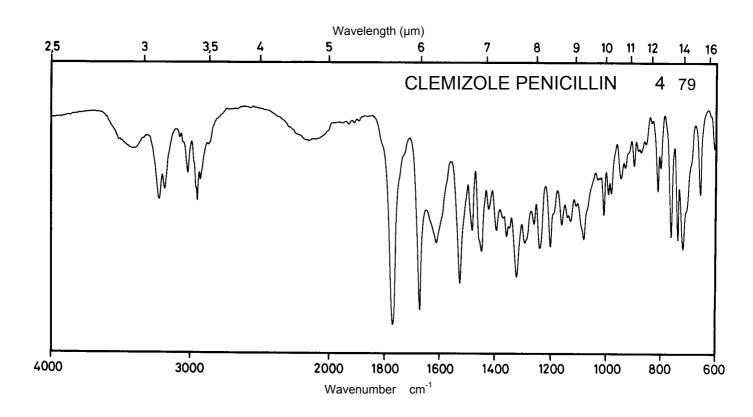




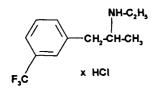


267.7


Concentration 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 263 nm			
E ^{1%} 1cm	20 21.9			
ε	546 590			

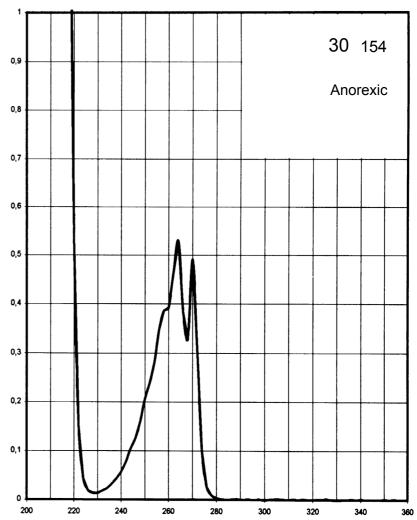

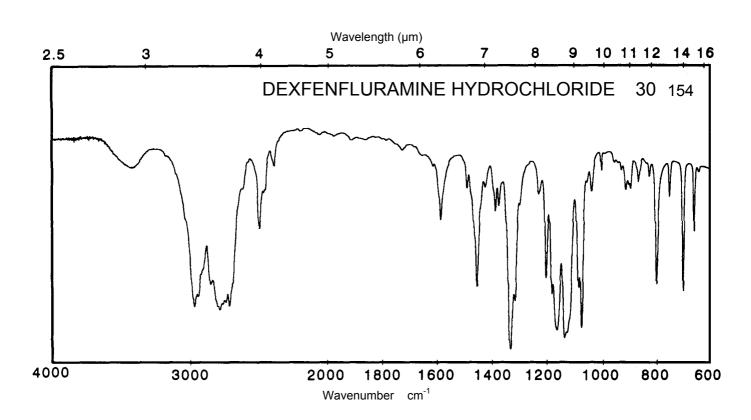




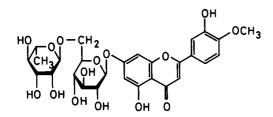


608.6


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			345 nm 266 nm 251 nm	372 nm 268 nm
E ^{1%} 1cm			336 292 292	229 440
3			20470 17760 17760	13900 26750

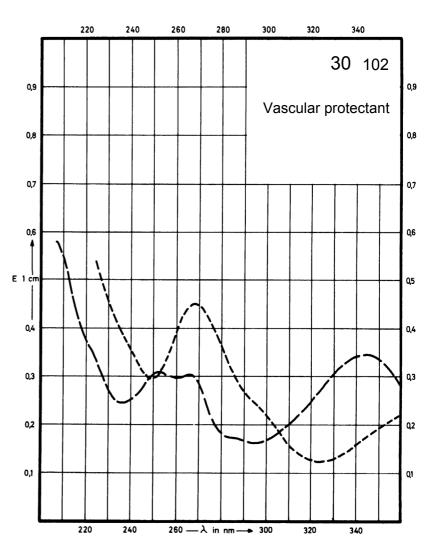

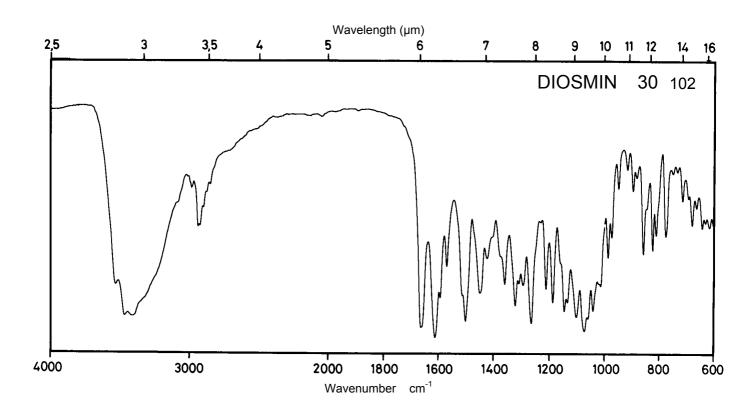




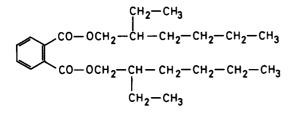


 M_r 390.6


Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm 224 nm			
E ^{1%} 1cm	32 218			
ε	1250 8520			

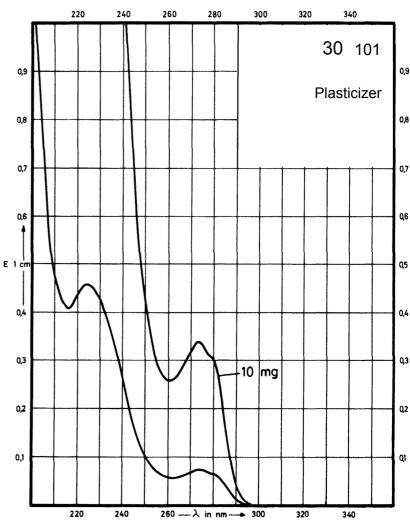

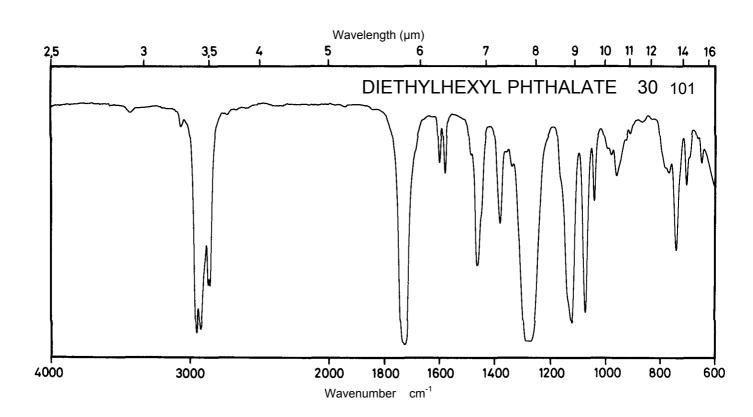



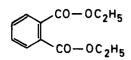



 M_r

222.2


Concentration 1.25 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm 224 nm			
E ^{1%} 1cm	56 375			
ε	1240 8340			

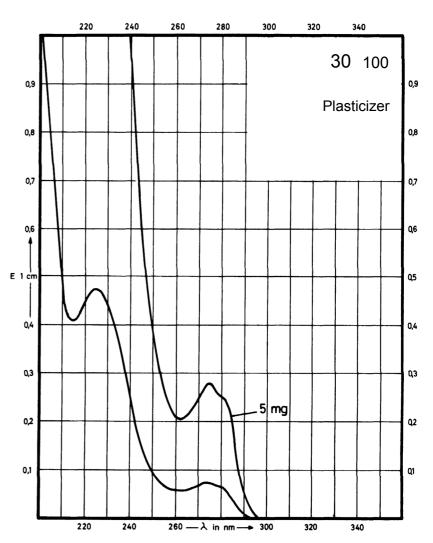

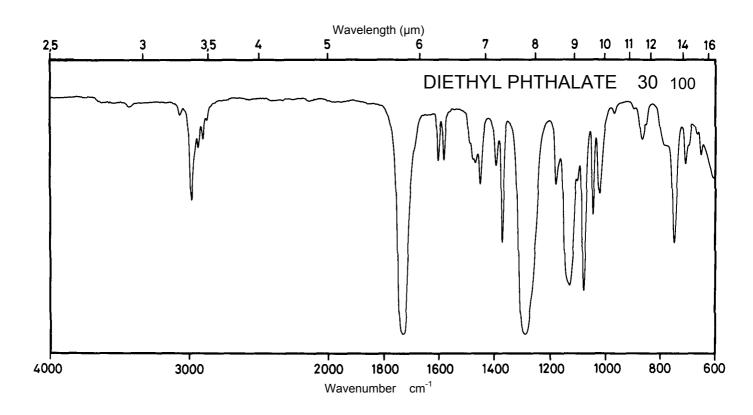



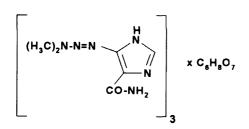



738.8

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	Decom- position observed	323 nm 224 nm	Decom- position observed
E ^{1%} _{1cm}			1116 436	
ε			82300 32200	

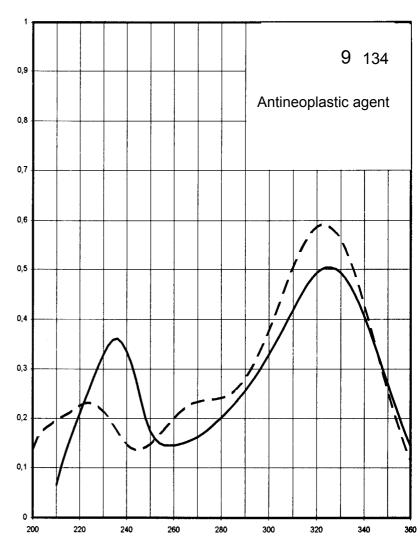

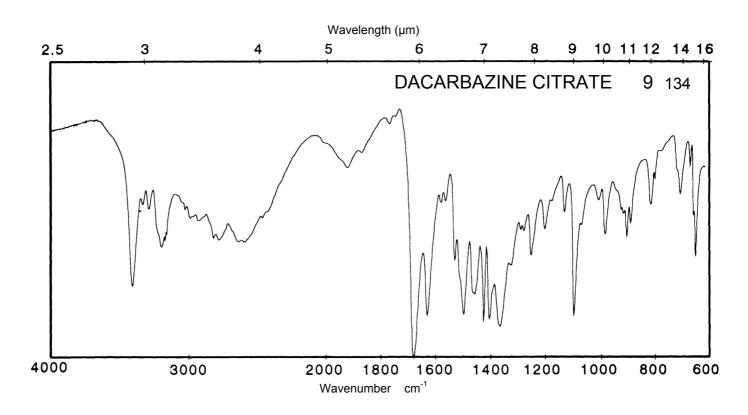




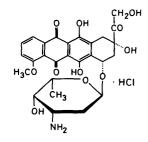


580.0


Concentration 1.05 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Syllibol				_
Maximum of	476 nm		478 nm	Decom-
	251 nm		253 nm	position
absorption	233 nm		232 nm	observed
10/_	215		205	
E ^{1%} 1cm	441		429	
1cm	660		634	
3	12500		11900	
	25600		24900	
	38800		36800	

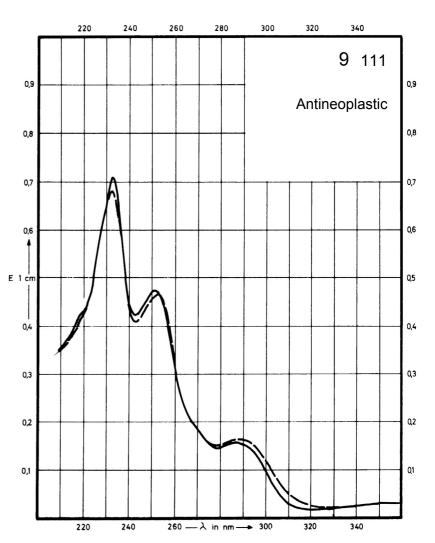

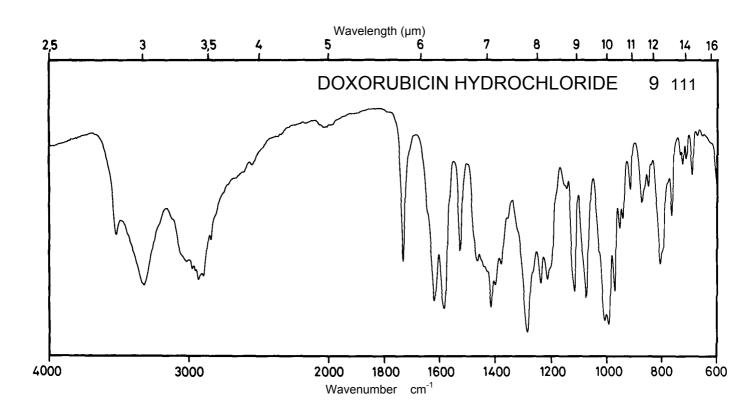

Name DAUNORUBICIN HYDROCHLORIDE

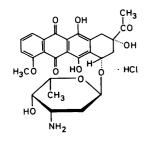



 M_r

564.0


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	494 nm 290 nm 233 nm		478 nm 290 nm 232 nm	590 nm 554 nm 252 nm
E ^{1%} 1cm	219 149 675		208 154 640	Decom- position observed
ε	12340 8380 38100		11700 8680 36080	

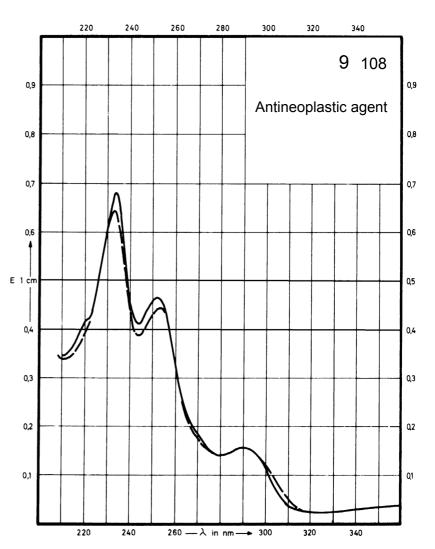

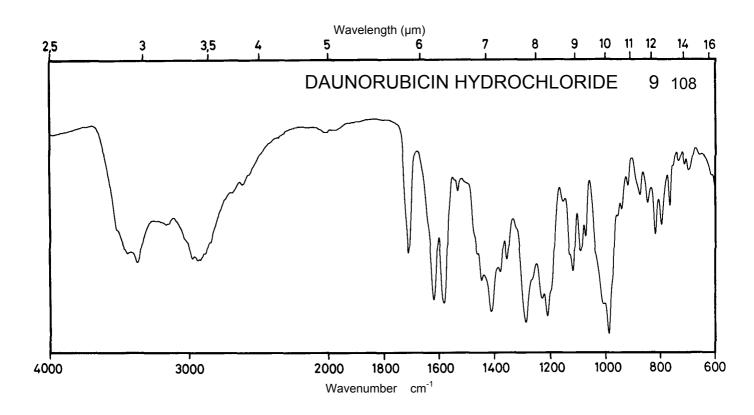




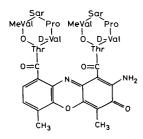


 M_r 1255.5


Concentration 2.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	440 nm 240 nm		442 nm 240 nm	343 nm 286 nm
E ^{1%} 1cm	187 260		174 251	Decom- position observed
ε	23500 32600		21850 31500	

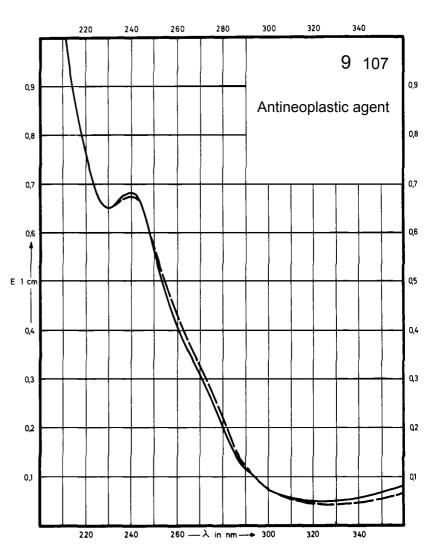

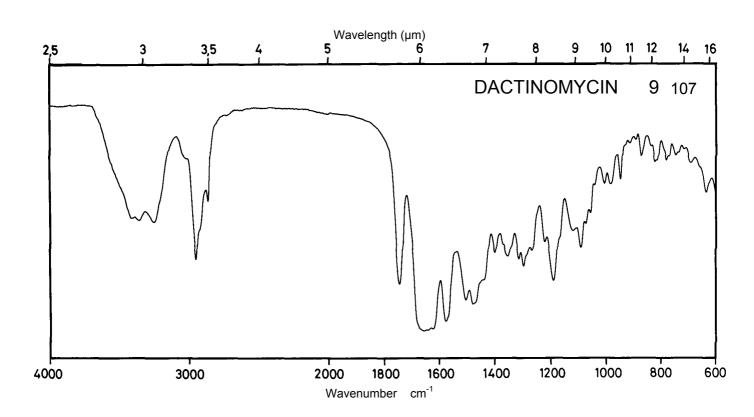




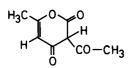


168.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	308 nm 224 nm		307 nm 224 nm	291 nm 230 nm
E ^{1%} 1cm	692 692		800 694	552 1150
ε	11630 11630		13460 11670	9290 19340

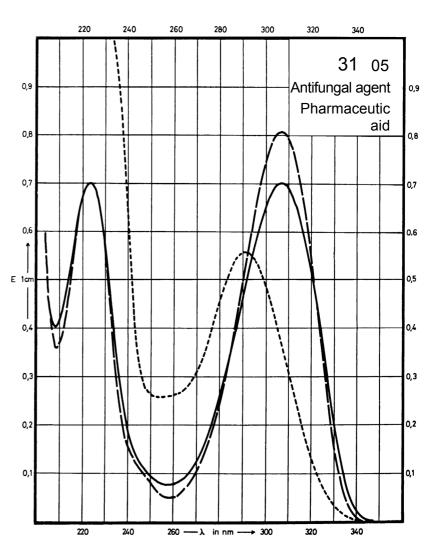

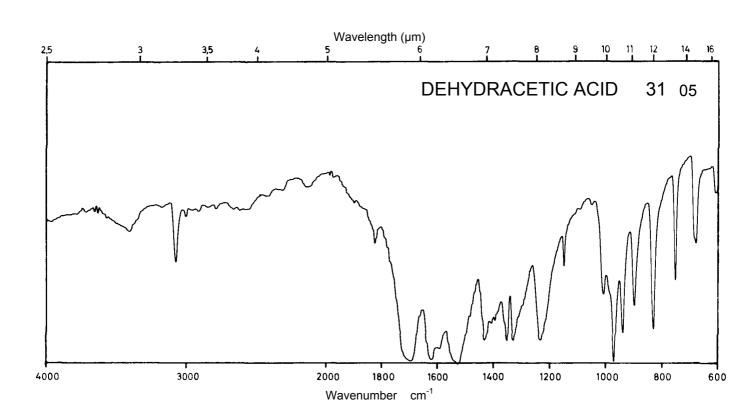




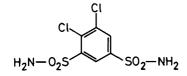


 M_{r} 305.2


Concentration 2 mg / 100 ml 6 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 293 nm 284 nm 295 nm 286 nm 293 nm 284 nm absorption E^{1%} 1cm 50 51 1520 1570 1750 1800 1100 1300 ε

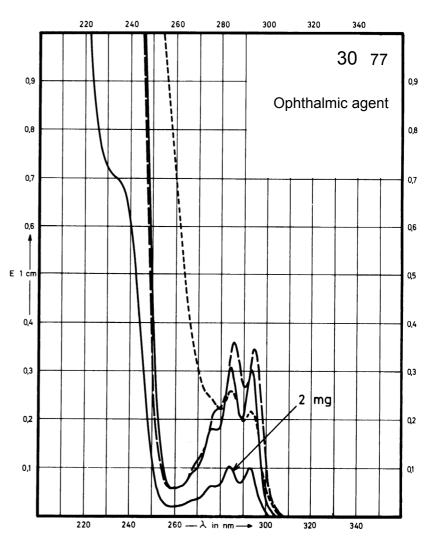

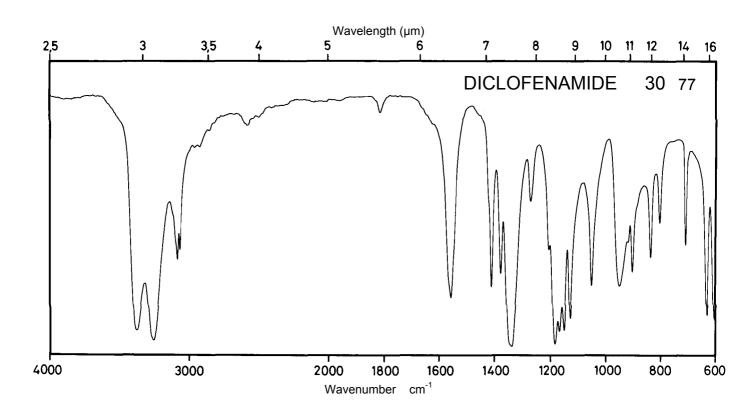



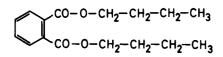



 M_r

278.4

Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm 224 nm			
E ^{1%} 1cm	45 303			
3	1250 8440			

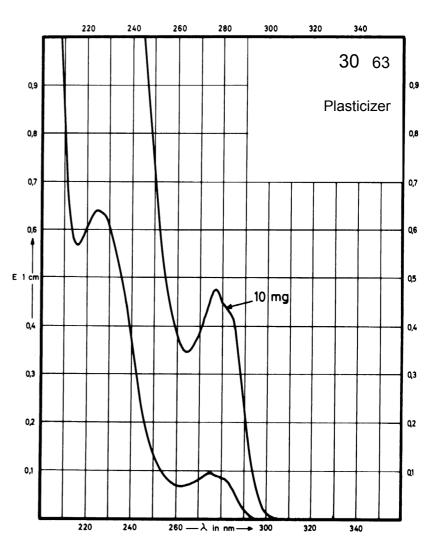

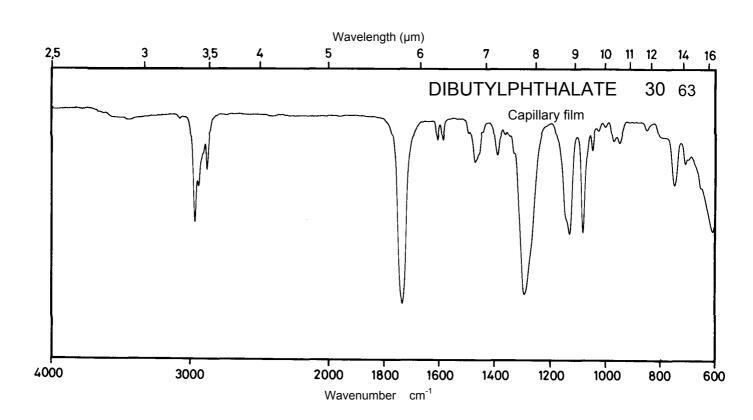



Name DISULFIRAM

 M_{r} 296.5


Concentration 1 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 325 nm 260 nm 245 nm 248 nm absorption E^{1%} 1cm 455 470 5820 25350 ε 13490 13930

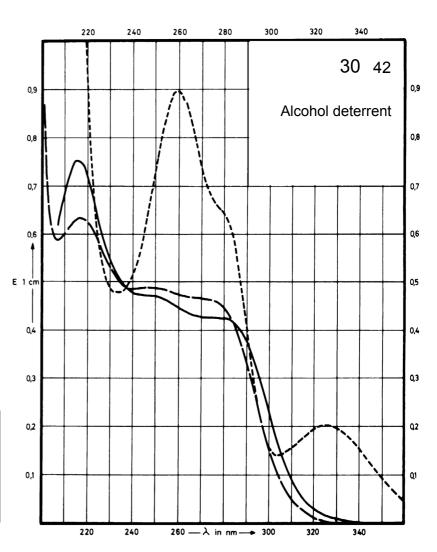

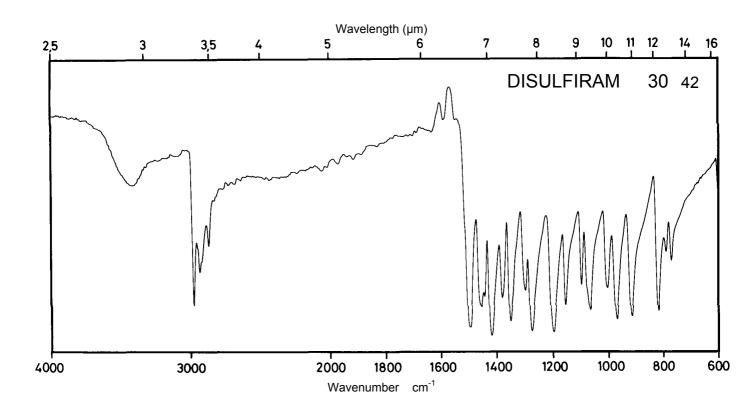



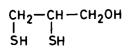



 M_r

Concentration 80 mg / 100 ml


124.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of absorption				
E ^{1%} 1cm				
ε				

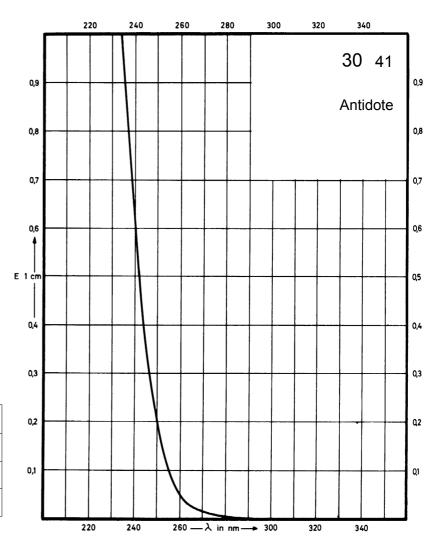

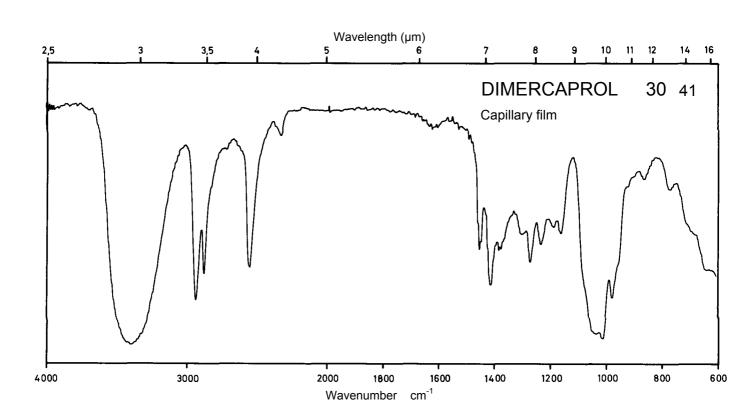





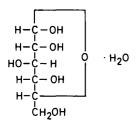


 M_r 198.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

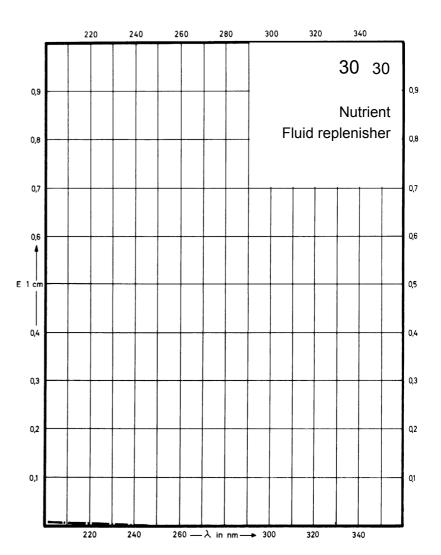

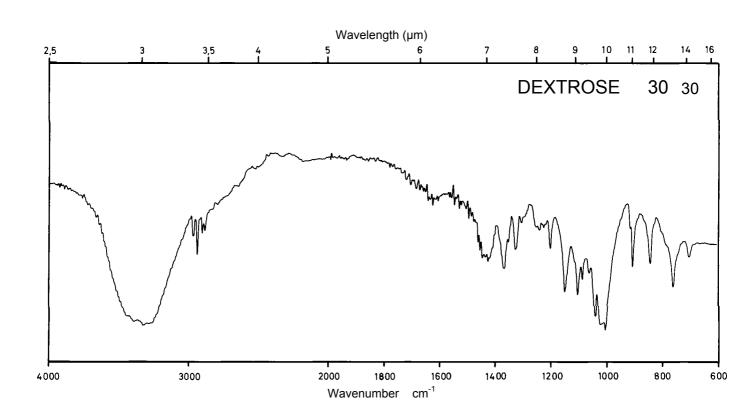





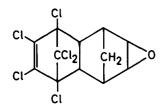

 M_r 380.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

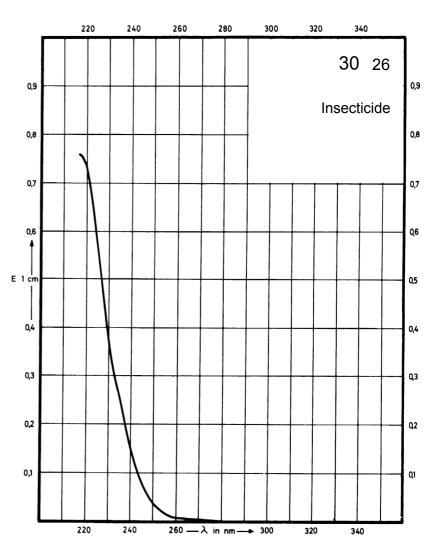

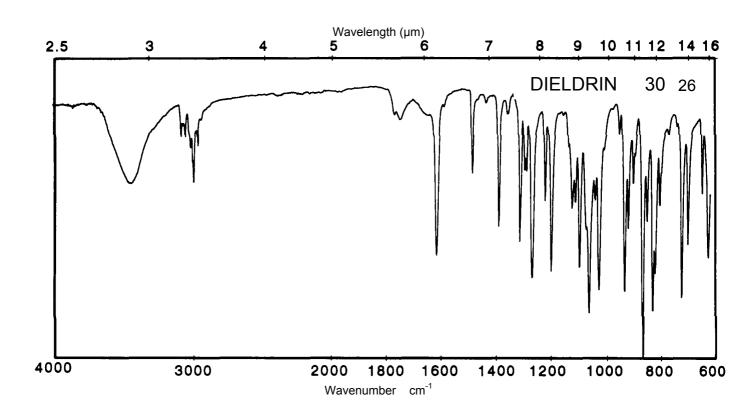





OH OH $N - (CH_2)_5 - NH - CO - (CH_2)_2 - CO - N - (CH_2)_5 - NH_2$ $CO - (CH_2)_2 - CO - NH - (CH_2)_5 - N - CO - CH_3$ OH


 M_{r} 560.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

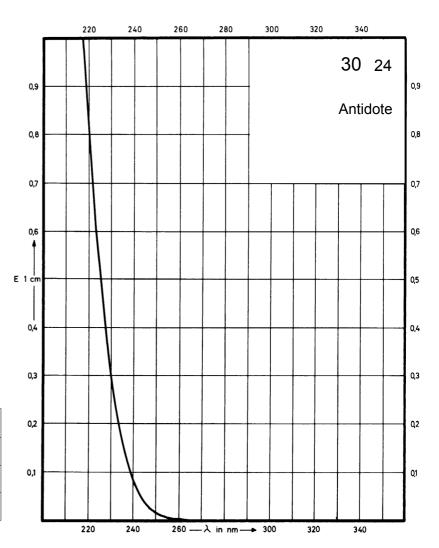

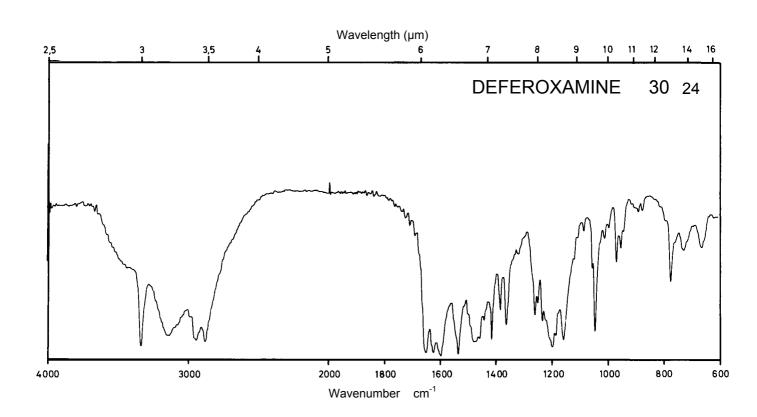


Name DESMETHYLDIA-ZEPAM

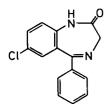


 M_r 270.7


Concentration 0.5 mg / 100 ml 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 228 nm		361 nm 283 nm 238 nm	342 nm
E ^{1%} 1cm	85 1300		155 500 1140	100
3	2300 35290		4190 13560 30810	2710

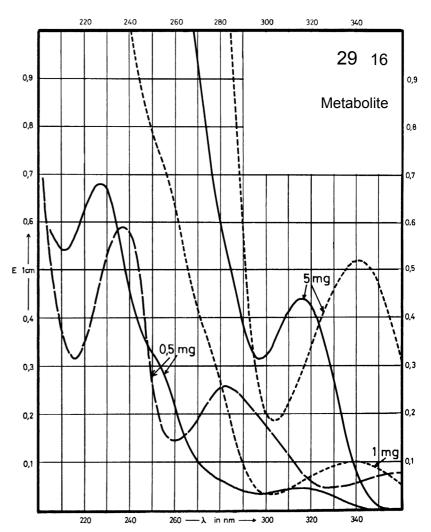

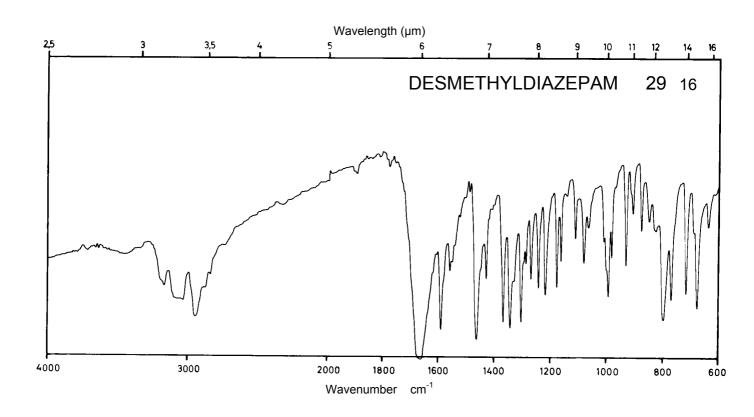


Name DOPAMINE HYDROCHLORIDE

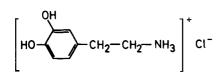



 M_{r}


189.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm		280 nm	
E ^{1%} 1cm	165		145	
3	3130		2750	

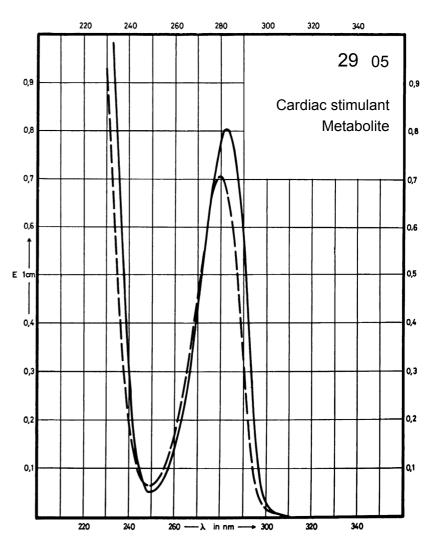


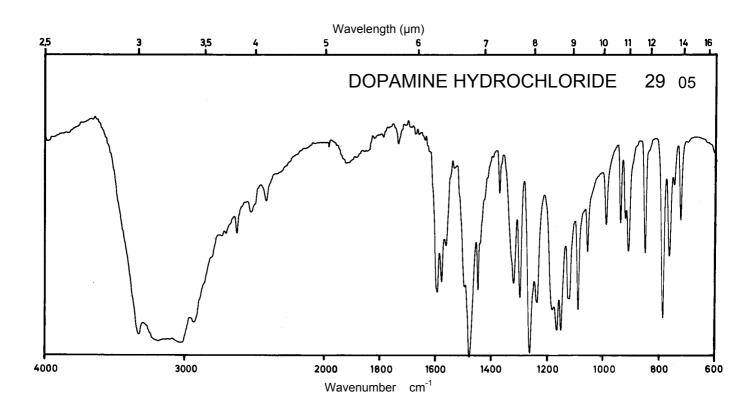

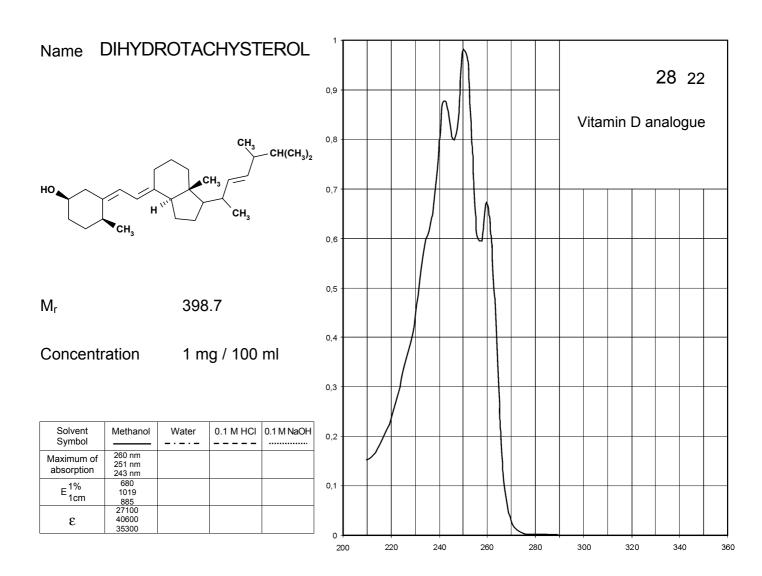

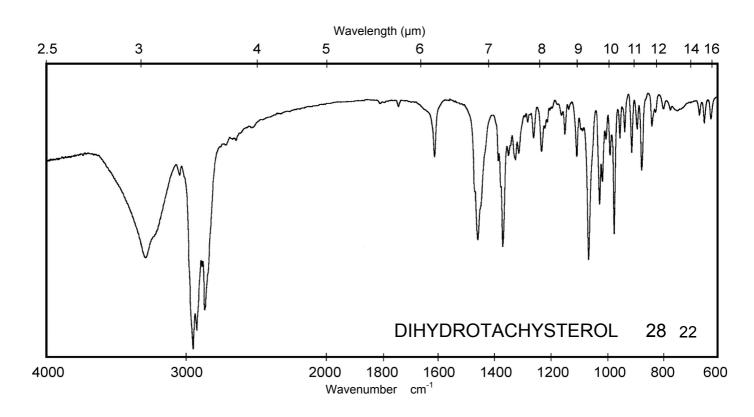




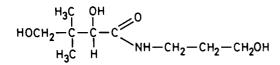


 M_r


Concentration 75 mg / 100 ml

205.3

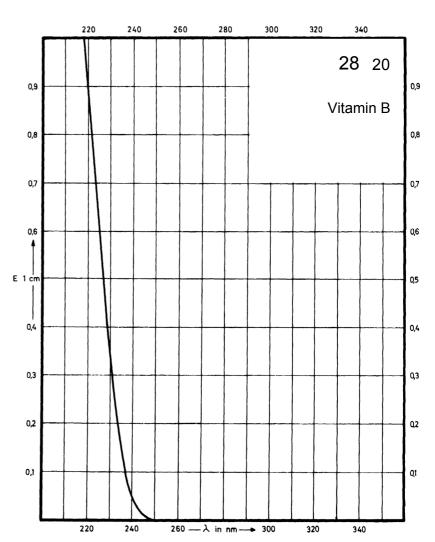

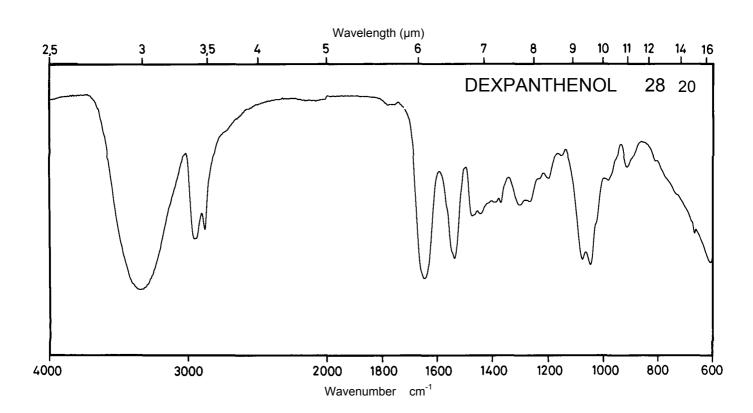



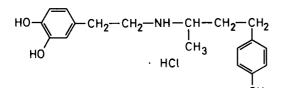


 M_{r}


337.8

Concentration


2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm 223 nm		278 nm 219 nm	Decom- position observed
E ^{1%} 1cm	142 405		122 385	
ε	4800 13700		4100 13000	

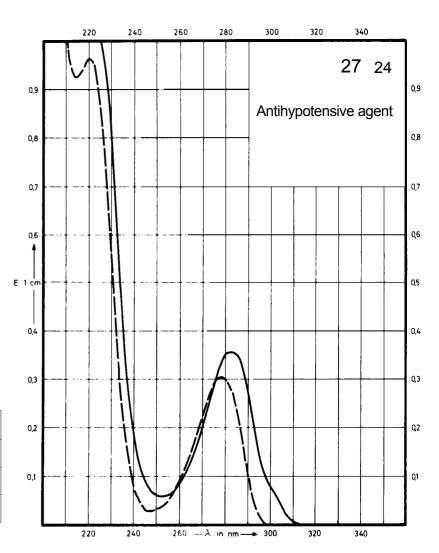

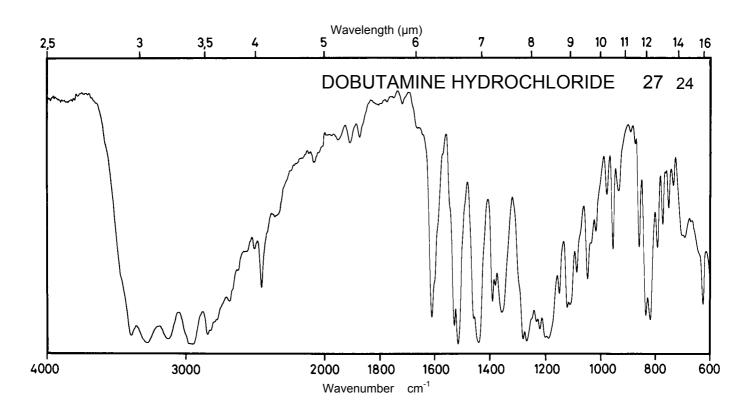





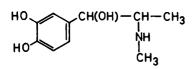

 M_{r} 197.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm		278 nm 220 nm	
E ^{1%} 1cm	141		122 270	
ε	2780		2410 5320	

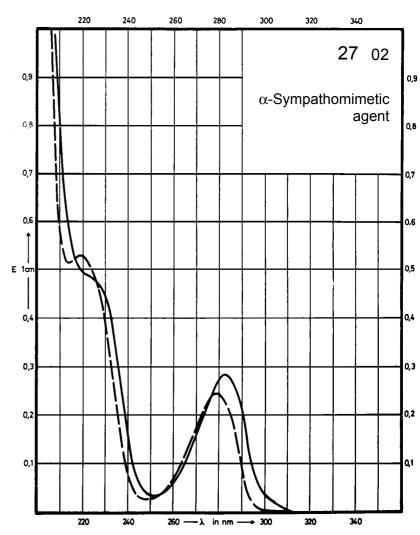

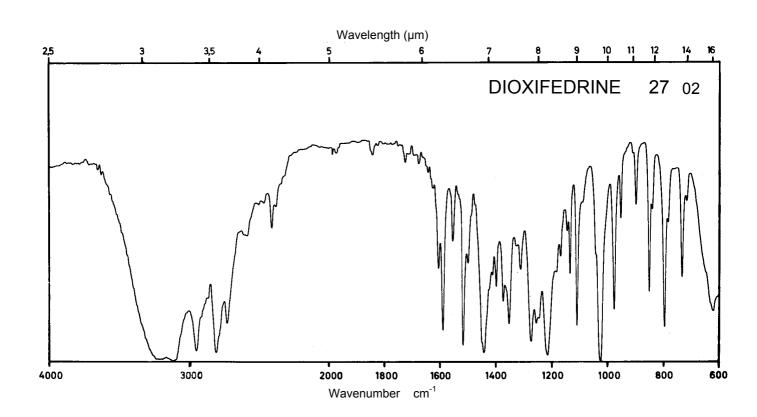


Name DOXYLAMINE SUCCINATE


$$\begin{array}{c|c} & \text{CH}_3 & \text{CH}_3 \\ & \text{C} - \text{O} - \text{CH}_2 - \text{CH}_2 - \text{N} \\ & \text{CH}_3 \\ & \text{C}_4 \text{H}_6 \text{O}_4 \\ \end{array}$$

 M_r 388.5

Concentration 2.5 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 261 nm 262 nm 261 nm absorption 109 232 109 ε 4240 9020 4240

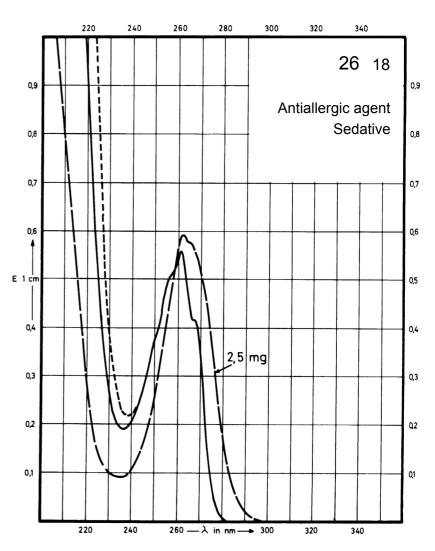

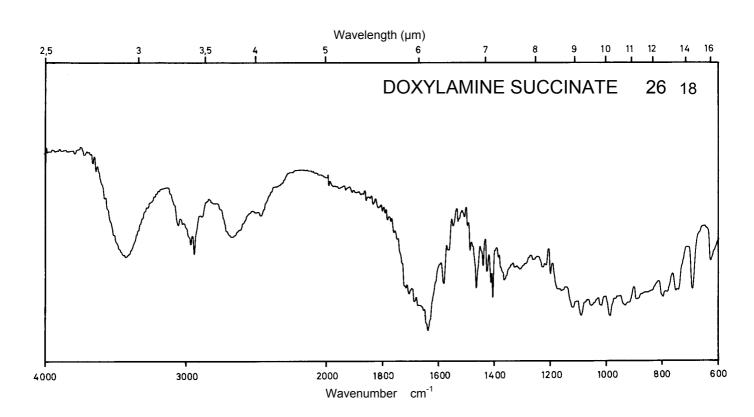


Name DIAZEPAM


CI-CH₃ 0 N-C C=N CH₂

 M_{r}


284.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm		360 nm 282 nm	310 nm
E ^{1%} 1cm	79		142 477	84
ε	2250		4040 13580	2390

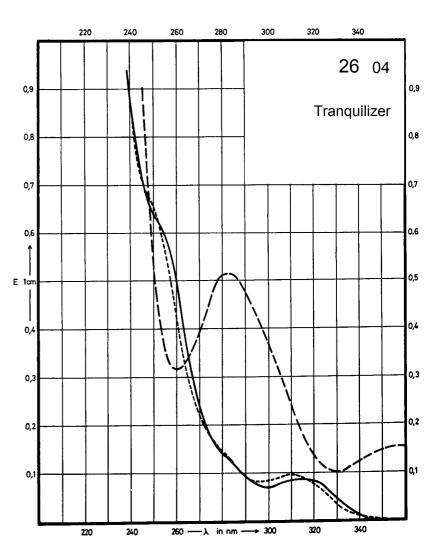

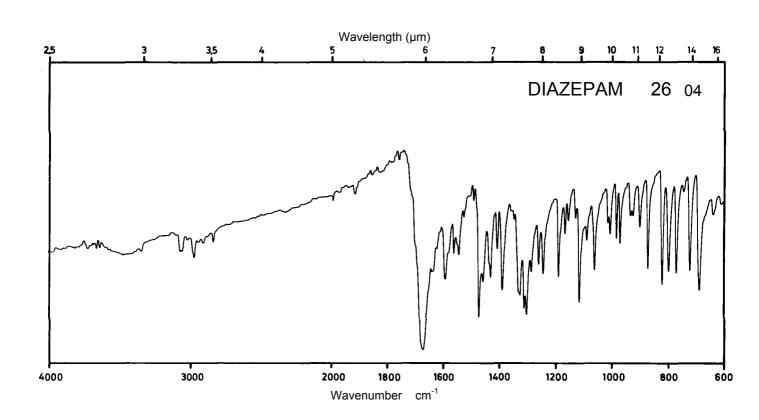




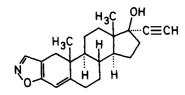


 M_r


337.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm			
E ^{1%} 1cm	338			
3	11400			

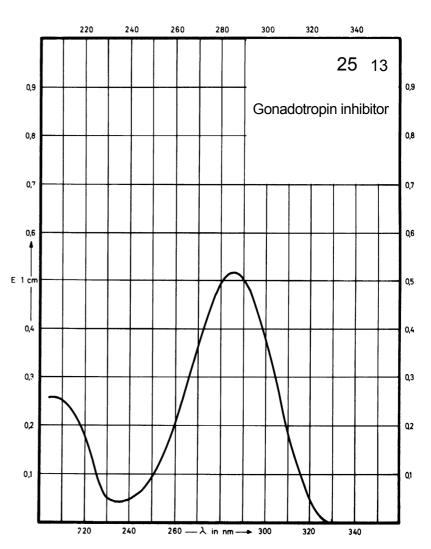

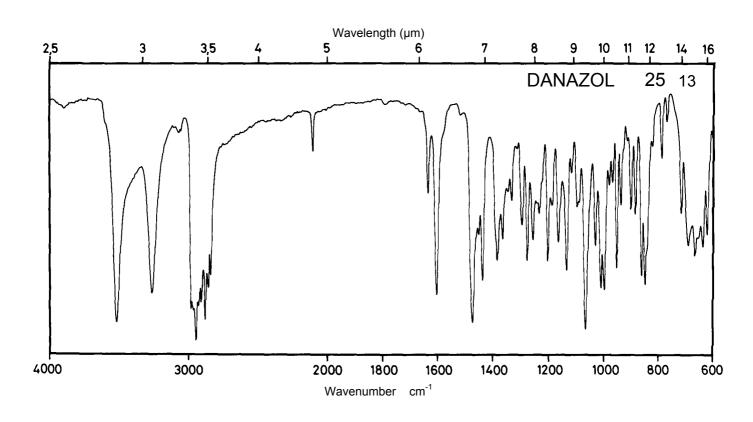





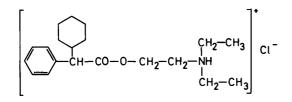


 M_{r} 353.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm		264 nm 258 nm 252 nm	
E ^{1%} 1cm	4.55 5.77 4.65		4.59 6.10 5.00	
3	160 205		162 215	

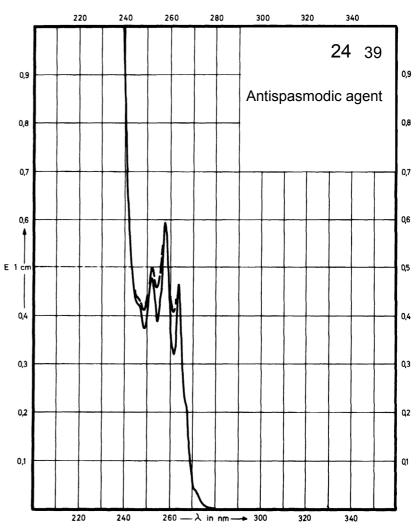

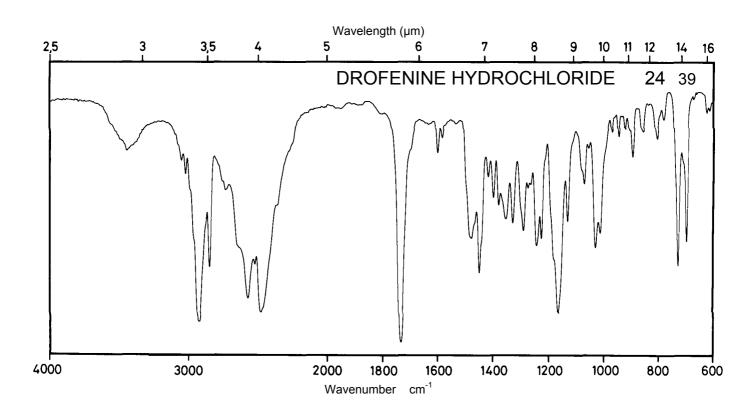




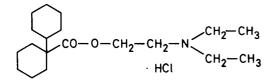


 M_r


346.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	216 nm		212 nm	
E ^{1%} 1cm	5.4		6.6	
ε	186		228	

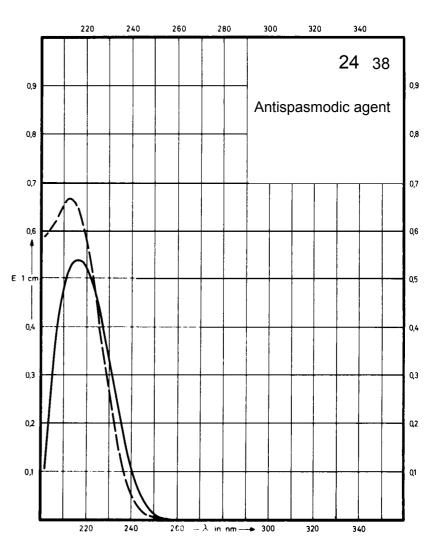

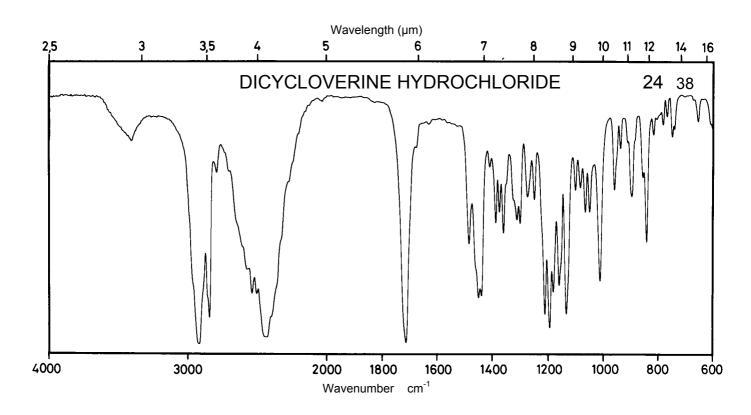

Name DANTROLENE SODIUM SALT

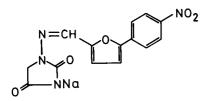


 M_{r}


336.2

Concentration


0.9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	381 nm 309 nm 226 nm		385 nm 312 nm 225 nm	397 nm 314 nm
E ^{1%} 1cm	838 395 368		774 407 402	748 487
ε	28170 13270 12370		26020 13690 13500	25150 16380

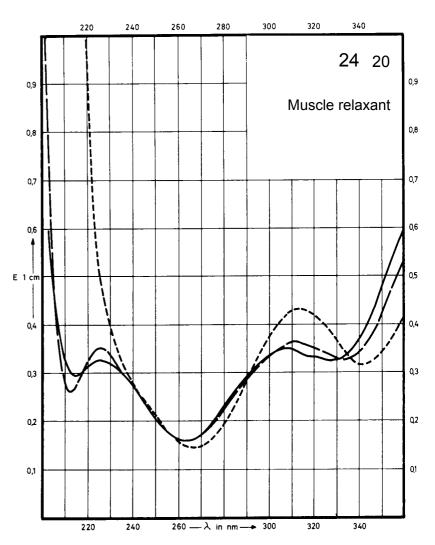

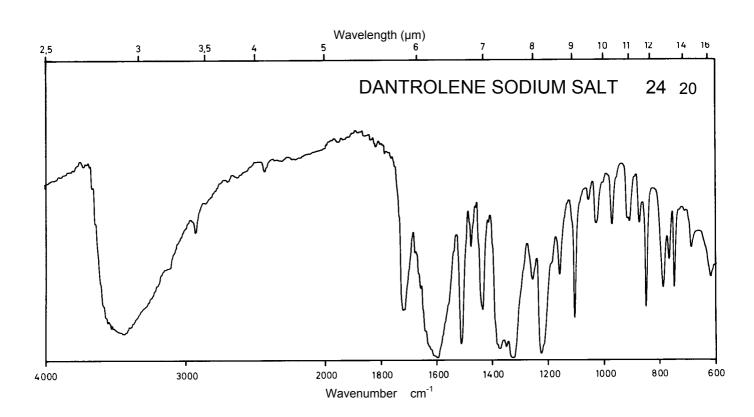


Name DIIODOTYROSINE

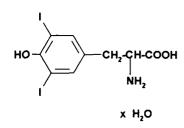


 M_{r}


451.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm 289 nm 221 nm	295 nm 220 nm	294 nm 221 nm	311 nm
E ^{1%} 1cm	65 67 660	58 623	55 625	125
3	2940 3030 29800	2620 28100	2500 28600	5650

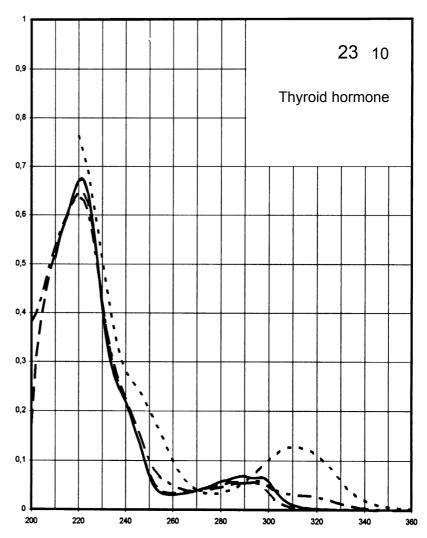

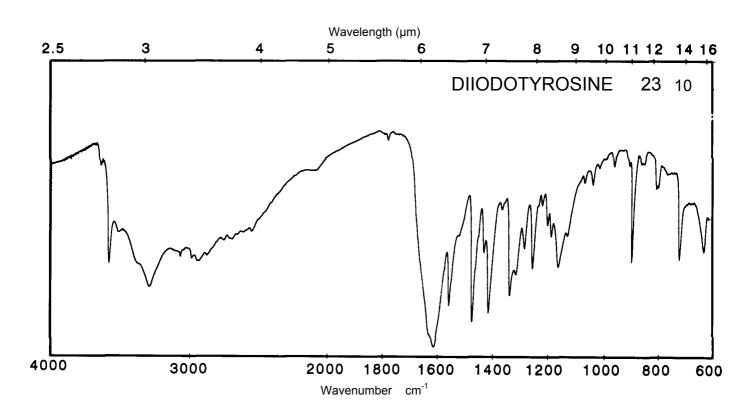





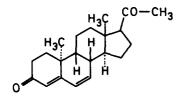


 M_{r} 312.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm			
E ^{1%} 1cm	840			
3	26250			

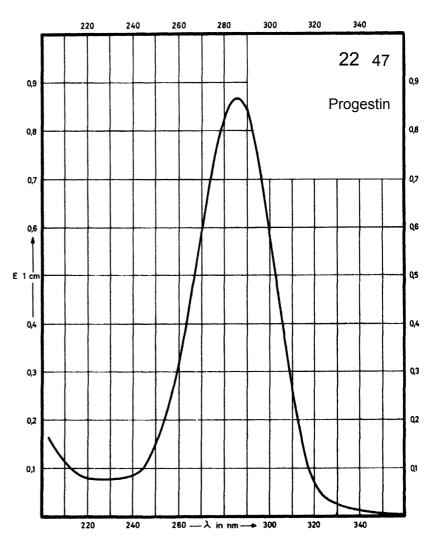

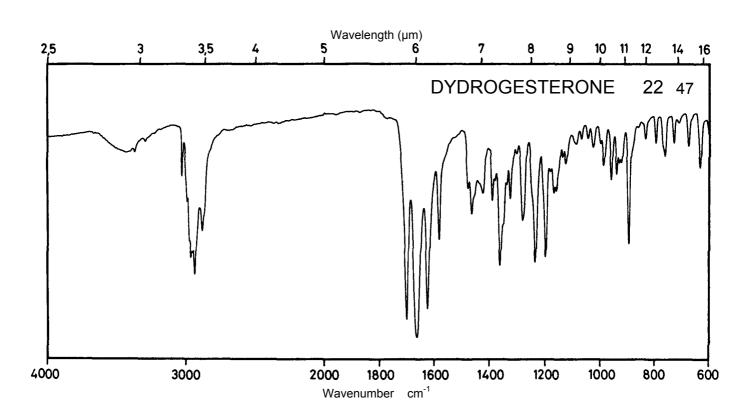




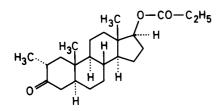


 M_{r}


360.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm			
E ^{1%} 1cm	0.7			
3	25			

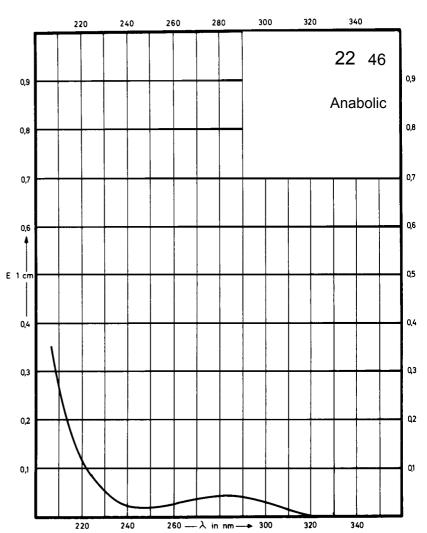

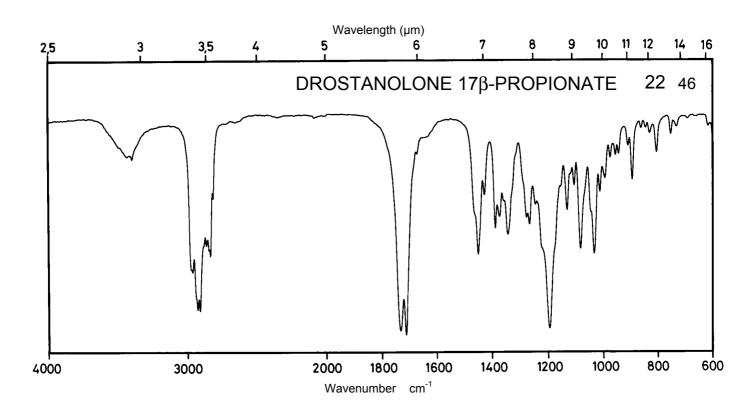





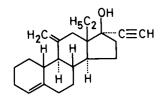


 M_r 286.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

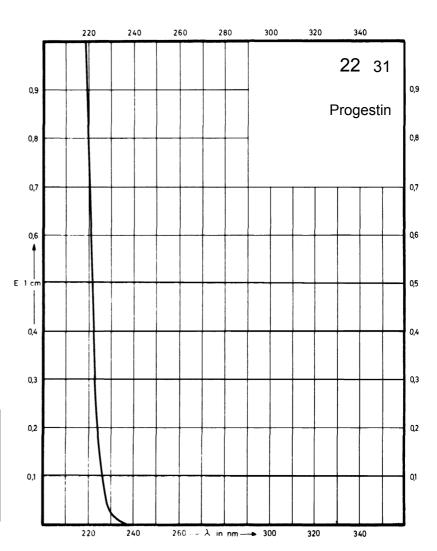

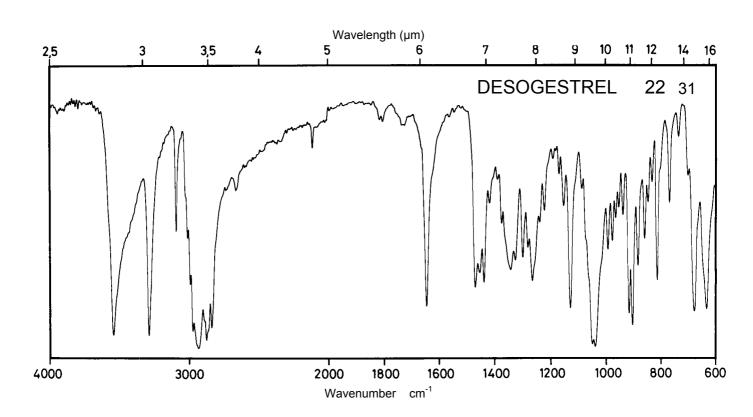


Name DIETHYLSTILBESTROL

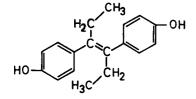


 M_{r}

268.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	236 nm			260 nm
E ^{1%} 1cm	409			452
ε	10980			12130

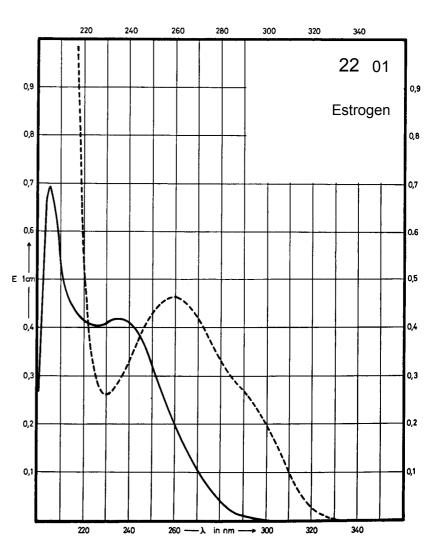

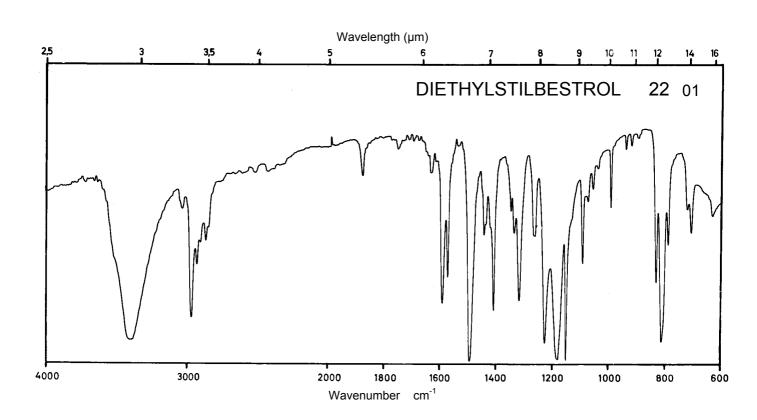



Name DOSULEPIN HYDROCHLORIDE

 M_{r}


331.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	305 nm 230 nm		302 nm 228 nm	
E ^{1%} 1cm	72 701		62 697	
3	2400 23250		2060 23150	

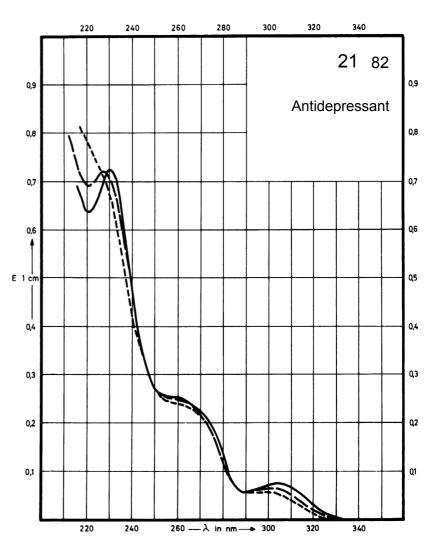

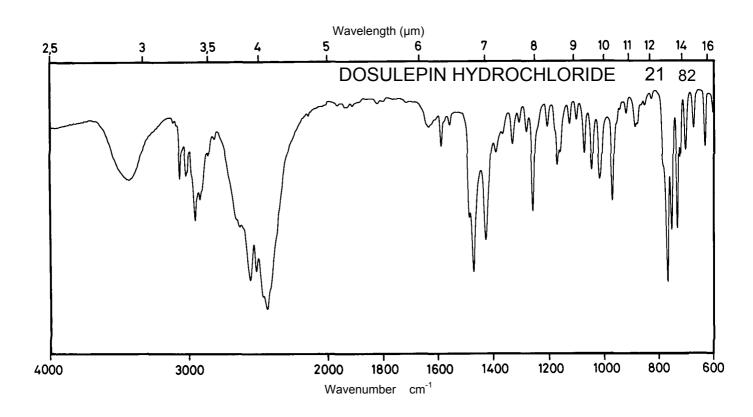



Name DEANOL 4-ACETAMIDO-BENZOATE

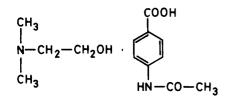

 M_r 268.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm		268 nm	261 nm
E ^{1%} 1cm	773		718	642
3	20740		19270	17240

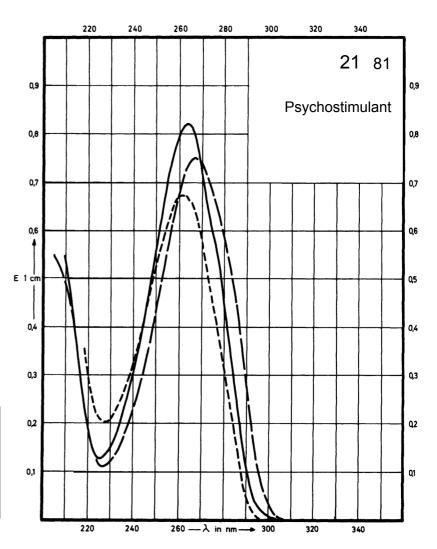

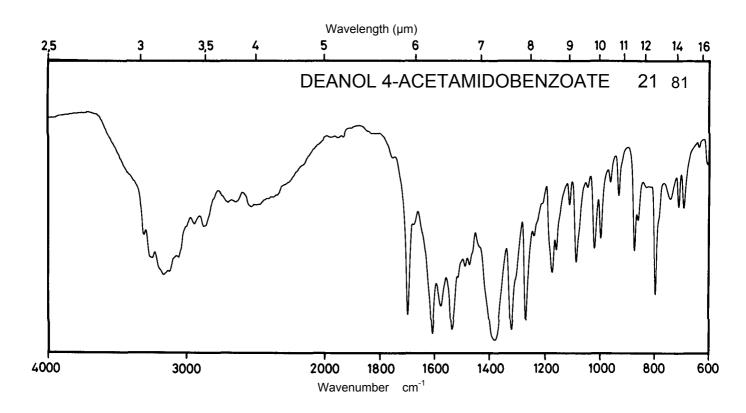



Name DIMETACRINE TARTRATE

$$\begin{array}{c|c} & & & & & \\ & & & & & \\ & & & & \\ & & & & \\ & & & & \\ & & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & \\ & & \\ & & \\ & & \\ & & \\ & & \\ & & \\ & & \\ & & \\ & & \\ & & \\ & \\ & & \\ &$$


 M_{r} 444.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	286 nm	286 nm	286 nm	
E ^{1%} 1cm	361	348	351	
3	16060	15480	15580	

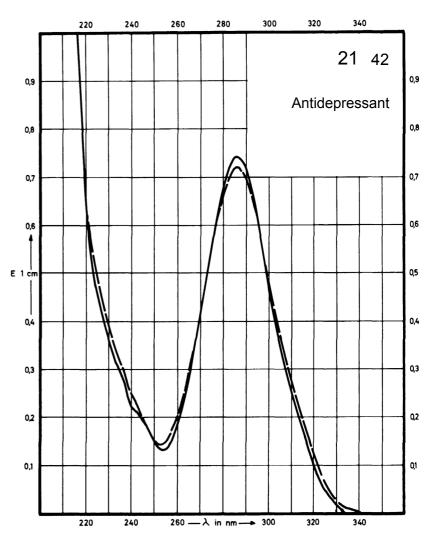

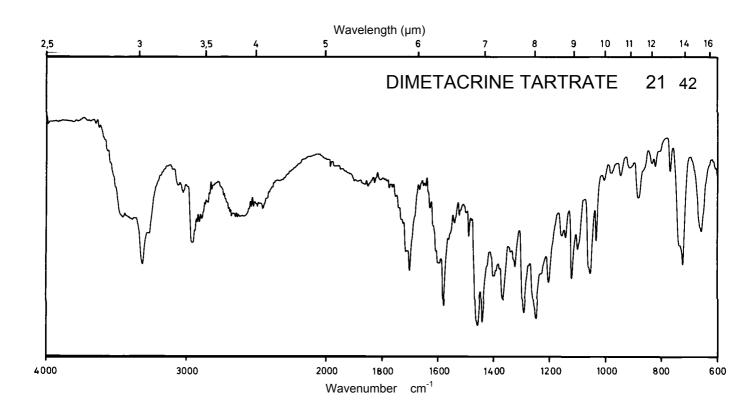

Name DROPERIDOL

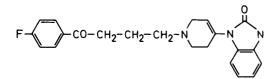



 M_{r}

379.4


Concentration 1 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm 233 nm		276 nm 247 nm	286 nm 246 nm
E ^{1%} 1cm	181 432		191 397	195 515
ε	6870 16370		7250 15060	7400 19540

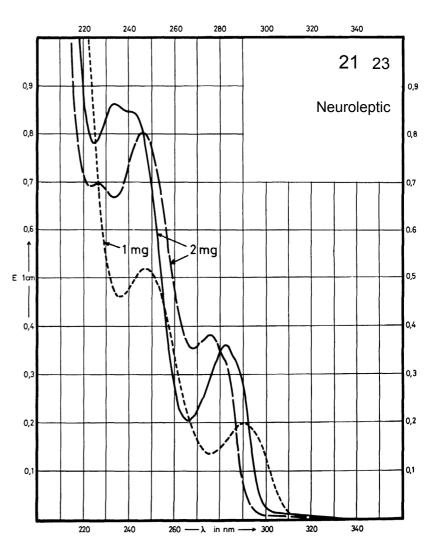

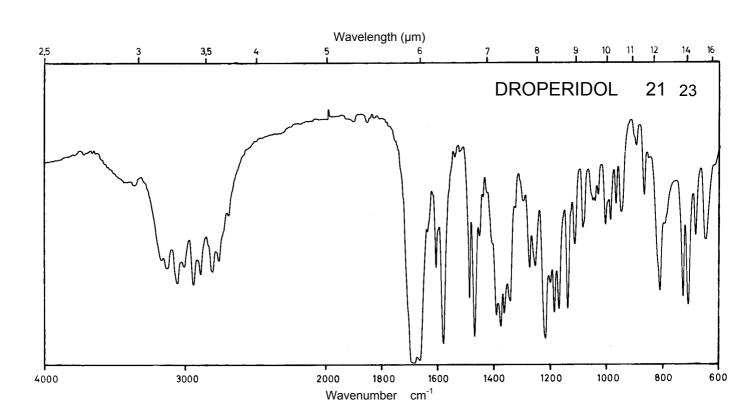


Name DIBENZEPIN HYDROCHLORIDE

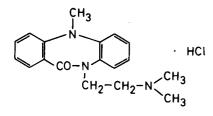


 M_{r}


331.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	225 nm		221 nm	
E ^{1%} 1cm	840		830	
3	27870		27540	

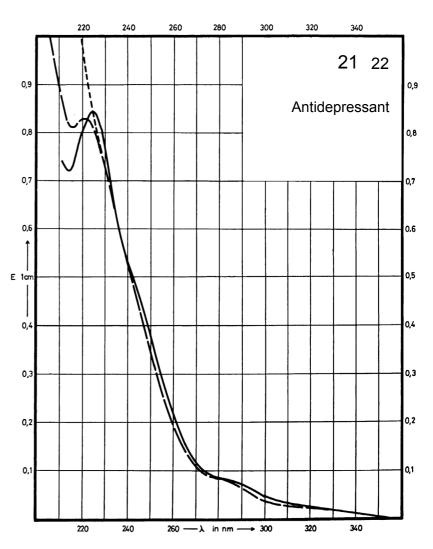

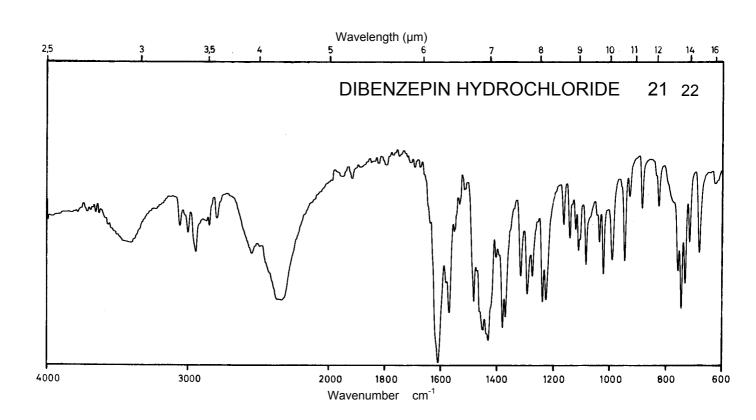


Name DESIPRAMINE HYDROCHLORIDE

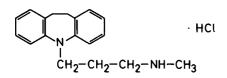


 M_{r}


302.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	250 nm		251 nm	254 nm
E ^{1%} 1cm	296		271	266
3	8950		8210	8060

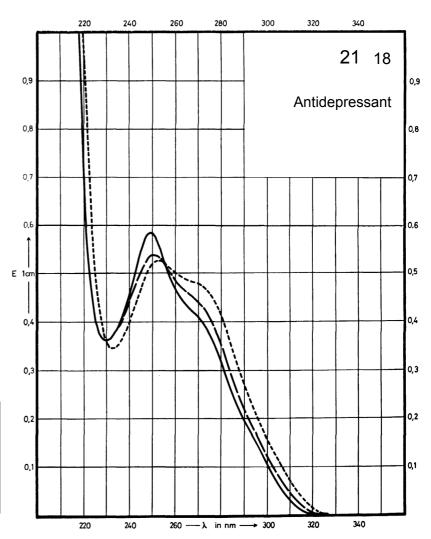

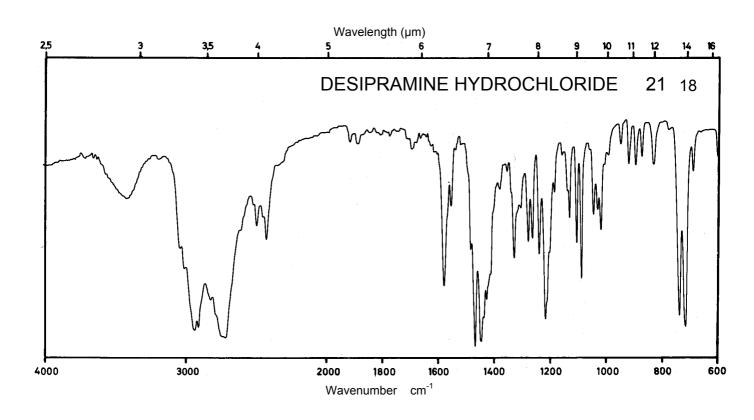


Name DOXEPIN HYDROCHLORIDE

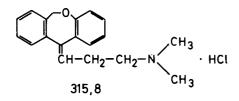


 M_{r}


315.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm		292 nm	
E ^{1%} 1cm	108		92	
ε	3410		2910	

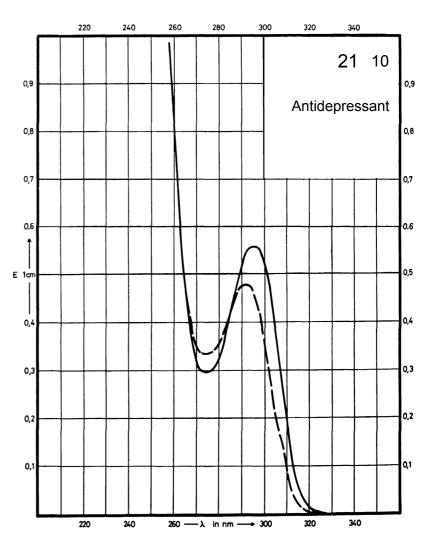

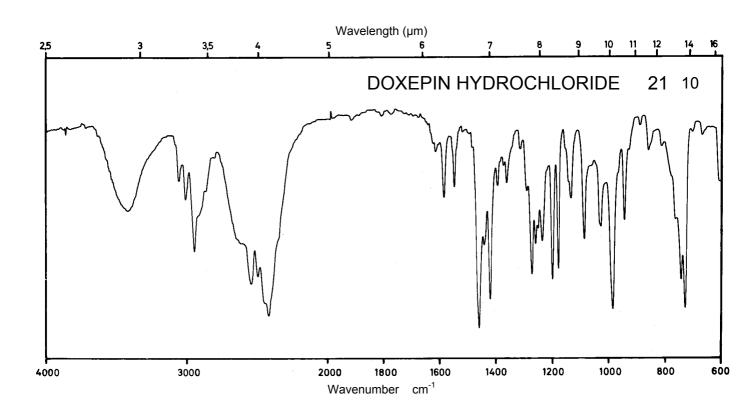





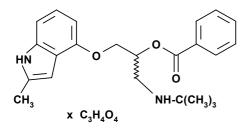


 M_r 484.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm 222 nm	267 nm 220 nm	267 nm 220 nm	267 nm
E ^{1%} 1cm	229 990	199 882	200 885	183
ε	11100 48000	9650 42800	9680 42900	8870

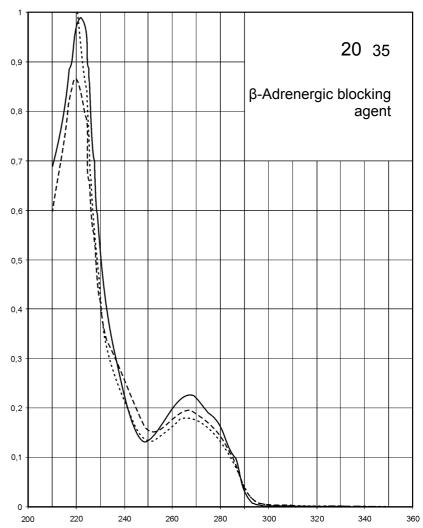

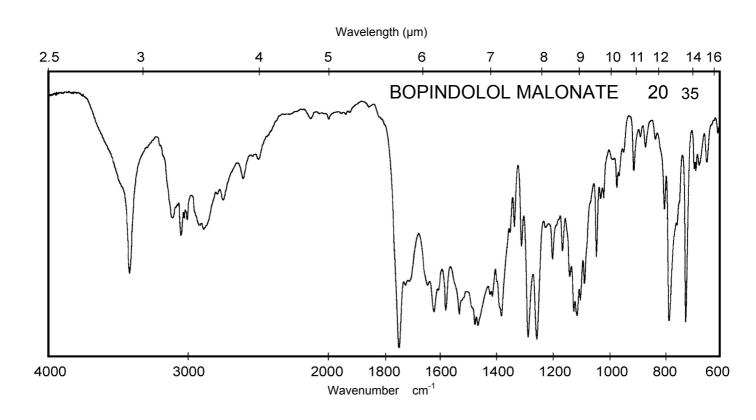




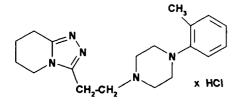


 M_{r}


361.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	235 nm	235 nm	234 nm	Decom- position observed
E ^{1%} 1cm	212	179	172	
3	7650	6450	6200	

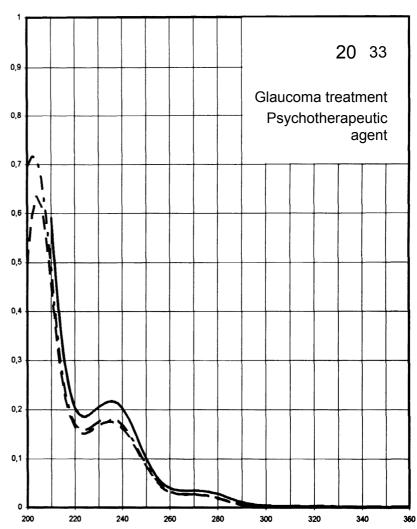

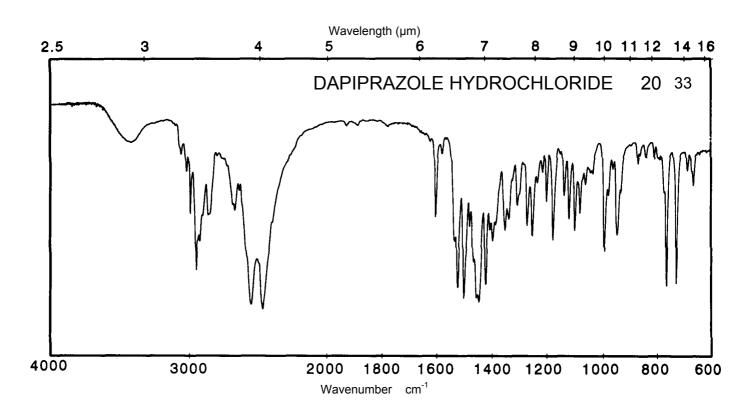

Name DIHYDROERGOT-AMINE MESILATE

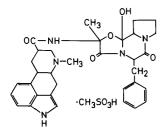


 M_{r}


679.8

Concentration


5.4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		279 nm	280 nm
E ^{1%} 1cm	100		94	91
ε	6800		6400	6180

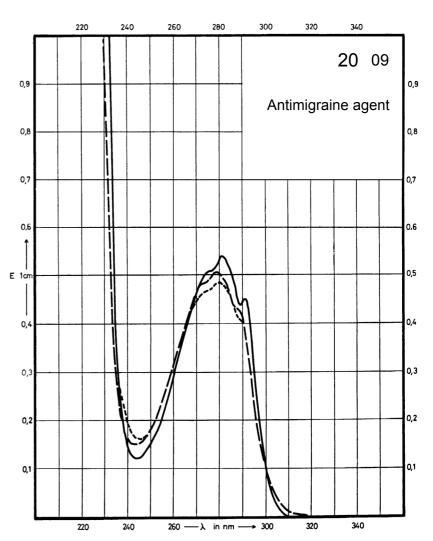

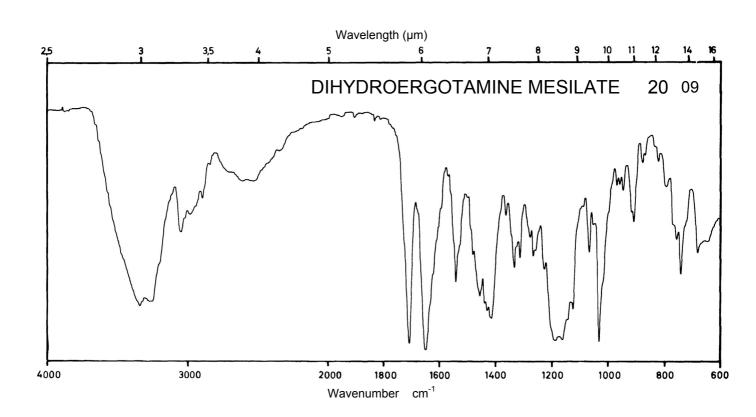



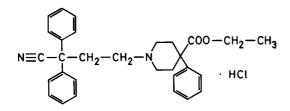



 M_{r}

489.1

Solvent Symbol	Methanol	Water	0.01 M HCI*	0.1 M NaOH
Maximum of absorption			264 nm 258 nm 252 nm	
E ^{1%} 1cm			10.2 12.9 10.6	
ε			501 632 520	

^{* 0.1} M HCl + propanol-2 (1+9)

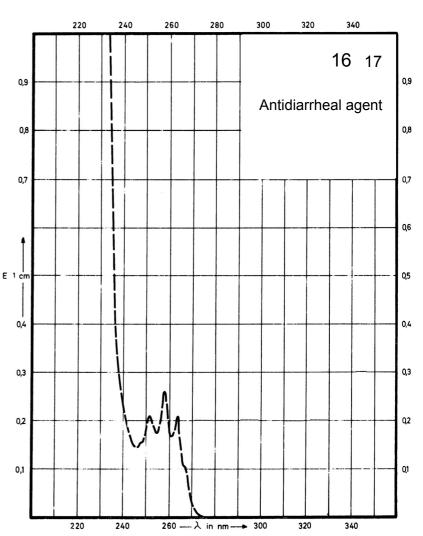

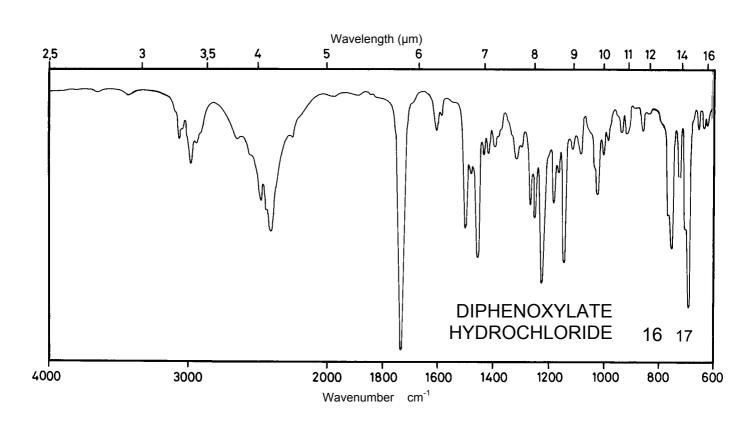





$$\begin{array}{c} \text{CH}_2-\text{CH}_3 \\ \text{CO}-\text{O}-\text{CH}_2-\text{CH}-\text{CH}_2-\text{CH}_2-\text{CH}_2-\text{CH}_3 \\ \text{CH}_2 \\ \text{CH}-\text{SO}_3\text{Na} \\ \text{CO}-\text{O}-\text{CH}_2-\text{CH}-\text{CH}_2-\text{CH}_2-\text{CH}_2-\text{CH}_3 \\ \text{CH}_2-\text{CH}_3 \end{array}$$

 M_r

444.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

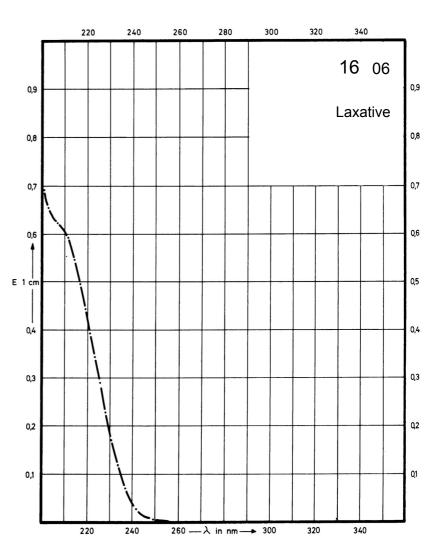

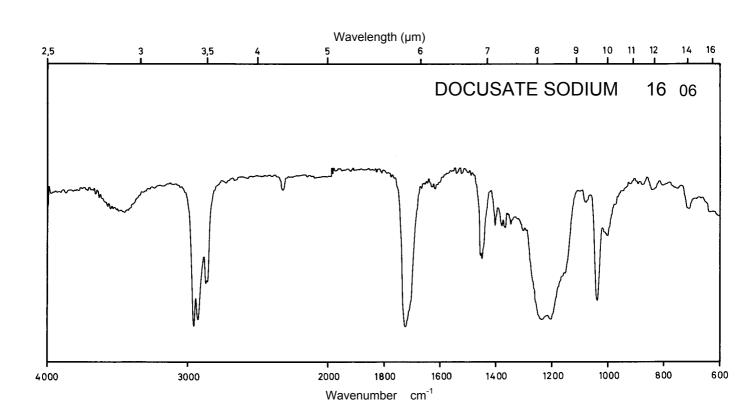



Name DEHYDROCHOLIC ACID

 M_r


402.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	286 nm			284 nm
E ^{1%} 1cm	5.0			6.3
ε	200			250

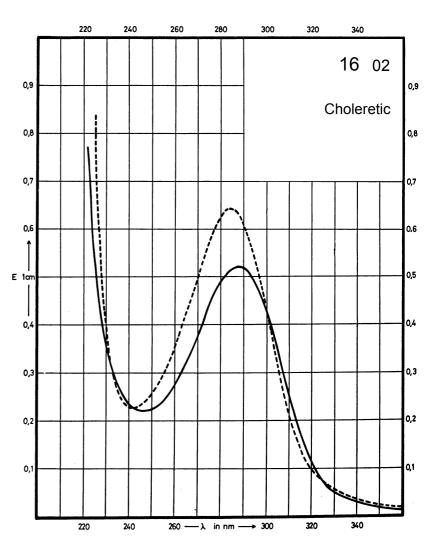

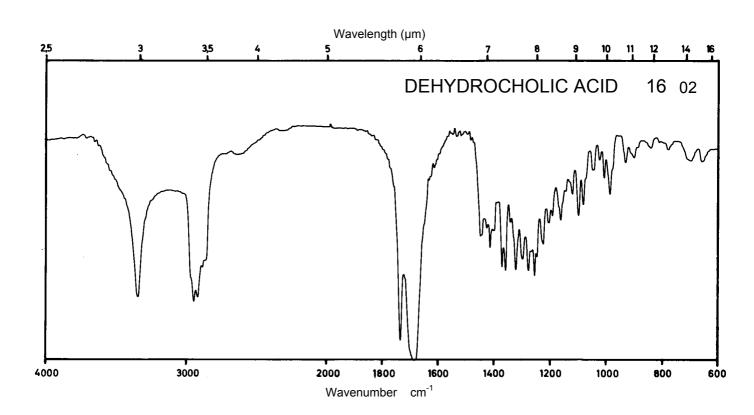




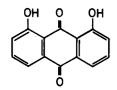


 M_r


240.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	500 nm 282 nm 250 nm			428 nm 278 nm
E ^{1%} 1cm	360 435 890			445 510
3	8650 10450 21380			10690 12250

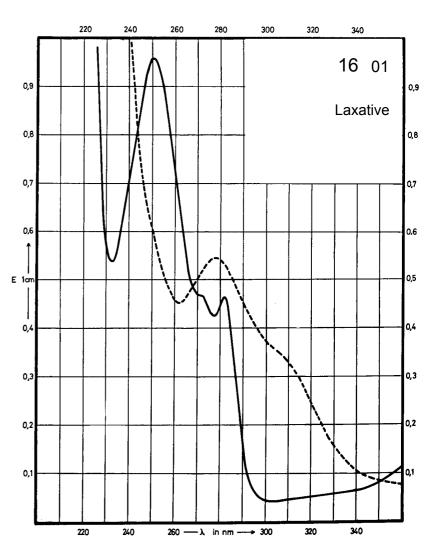

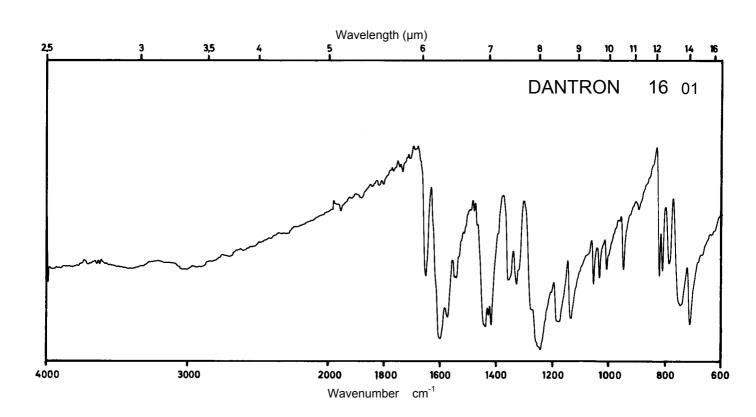




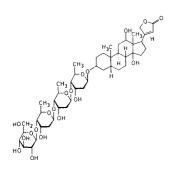


 M_{r}


943.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	218 nm			
E ^{1%} 1cm	153			
ε	14400			

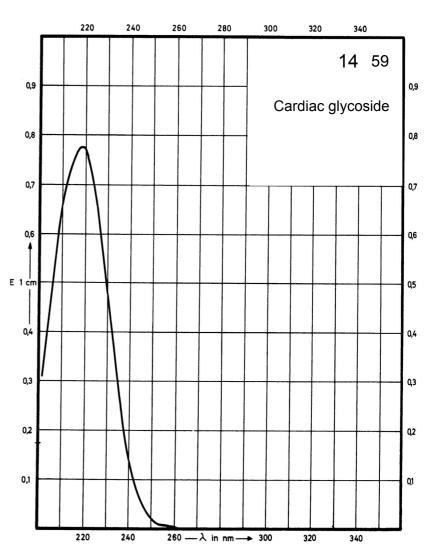

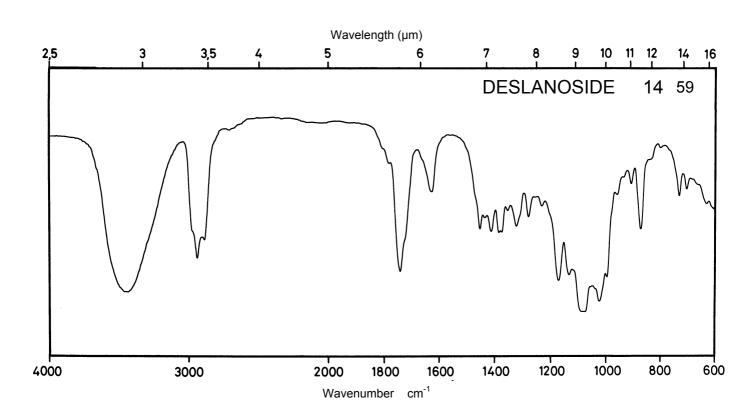




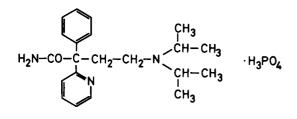


 M_{r}


437.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm		262 nm 267 nm	261 nm
E ^{1%} 1cm	182		300 299	187
ε	7980		13130 13080	8180

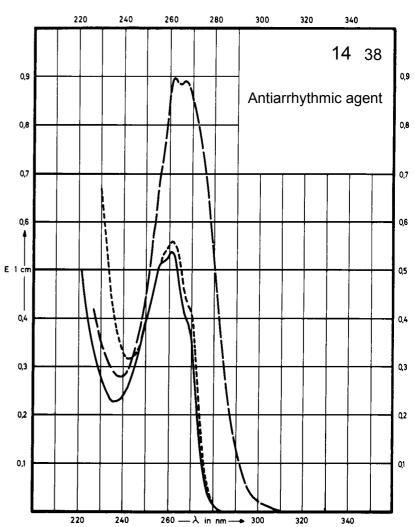

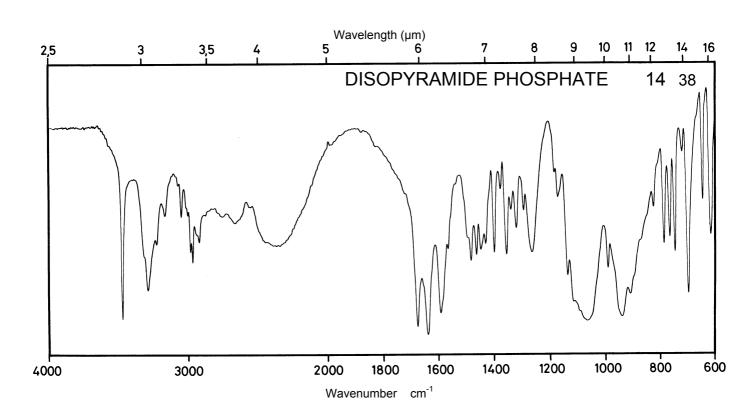


Name DILTIAZEM HYDROCHLORIDE

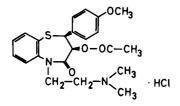


 M_{r}


451.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm	236 nm	236 nm	237 nm
E ^{1%} 1cm	537	522	533	555
ε	24200	23540	24040	25020

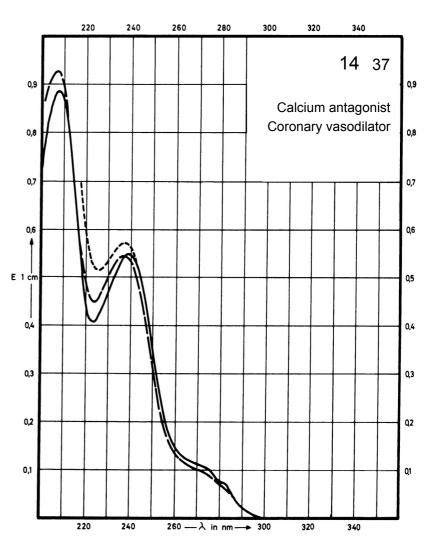

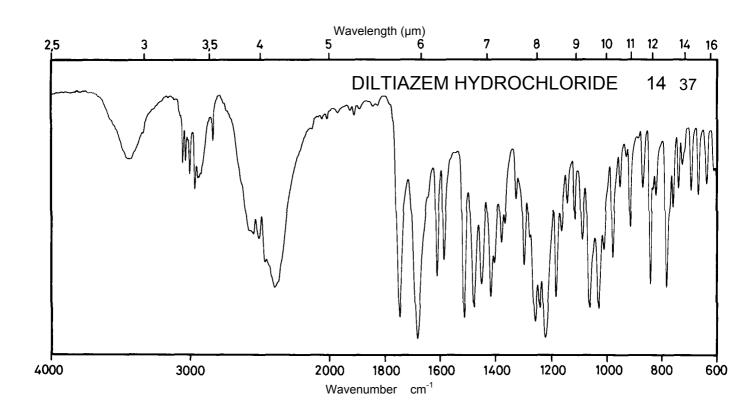





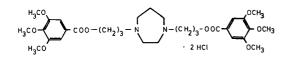


677.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm		265 nm	Decom- position observed
E ^{1%} 1cm	297		296	
3	20120		20060	

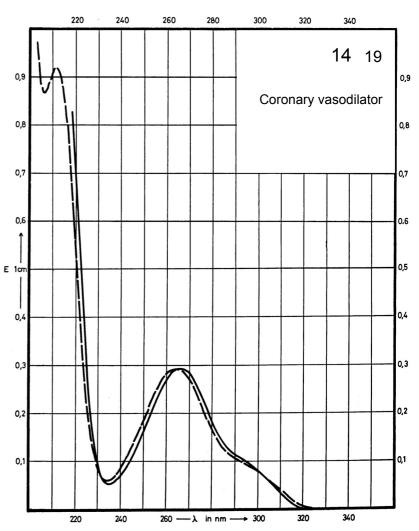

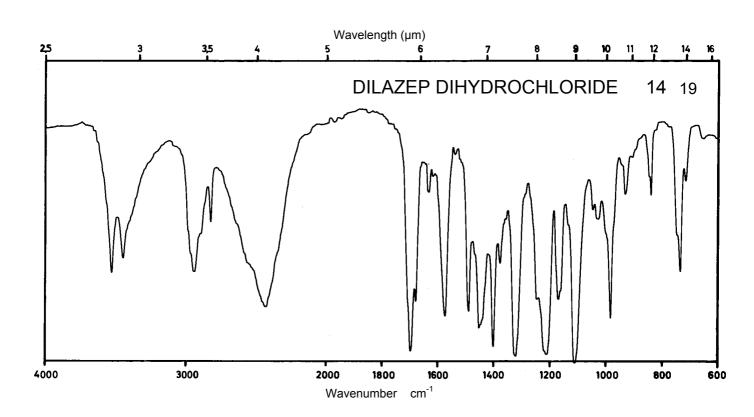


Name DIPYRIDAMOLE

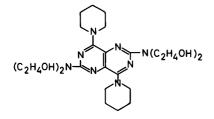


 M_{r}

504.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	408 nm 292 nm 228 nm		400 nm 283 nm 237 nm	414 nm 295 nm
E ^{1%} 1cm	165 570 520		134 546 554	155 558
3	8330 28760 26190		6760 27550 27950	7820 28160

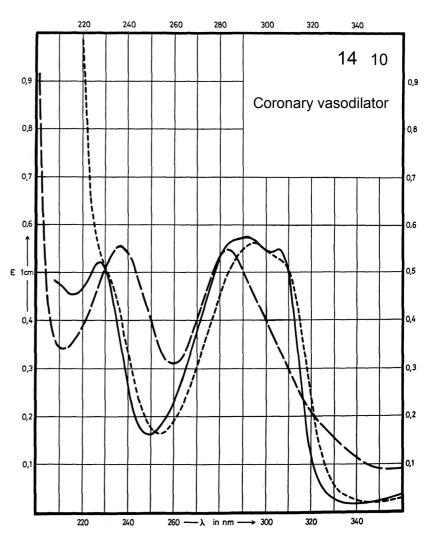

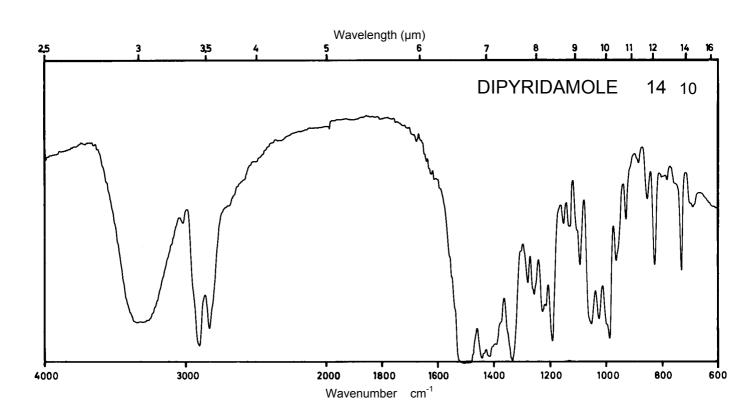


Name DIPROPHYLLINE


 $\begin{array}{c|c} \mathsf{CH}_3 & \mathsf{CH}_2 - \mathsf{CHOH} - \mathsf{CH}_2 \mathsf{OH} \\ \mathsf{CH}_3 & \mathsf{N} & \mathsf{N} \\ \mathsf{CH}_3 & \mathsf{CH}_3 \end{array}$

 M_{r}


254.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		272 nm	272 nm
E ^{1%} 1cm	344		365	360
ε	8740		9280	9150

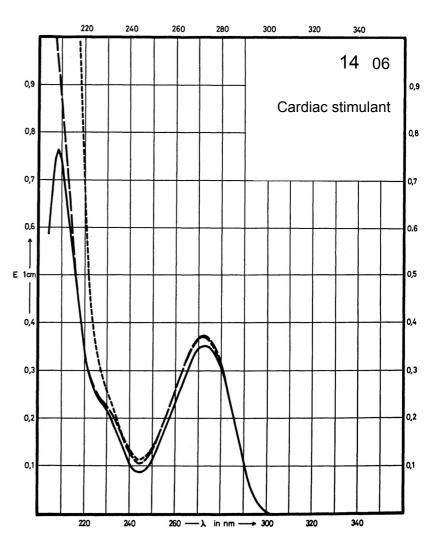

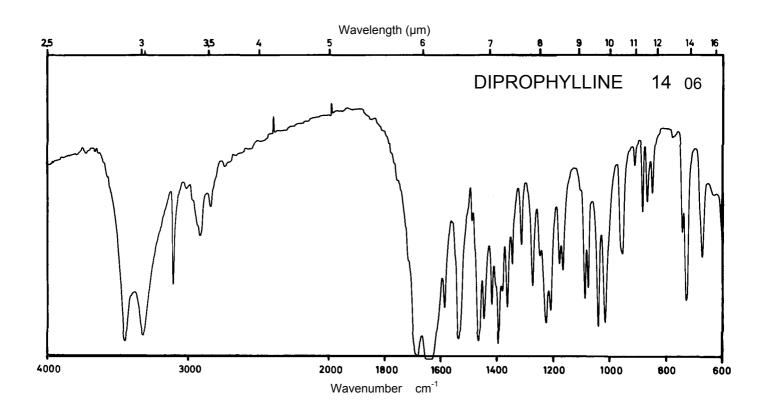



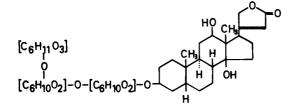



781.0


Concentration


2.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	218 nm		220 nm	
E ^{1%} 1cm	226		194	
ε	17650		15150	

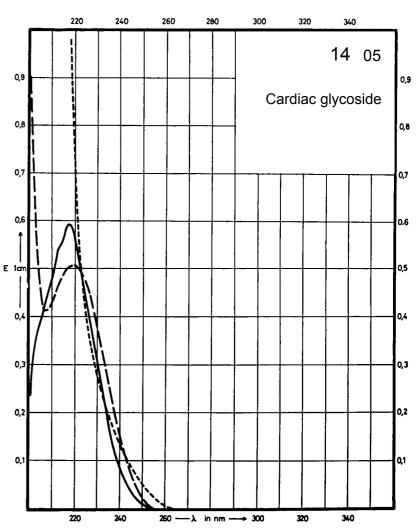


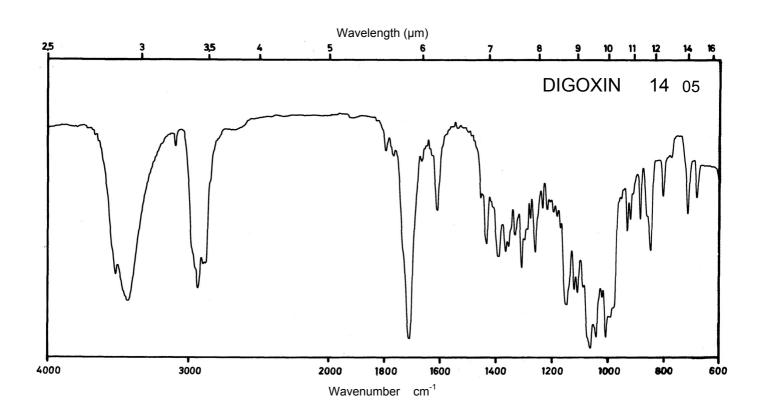

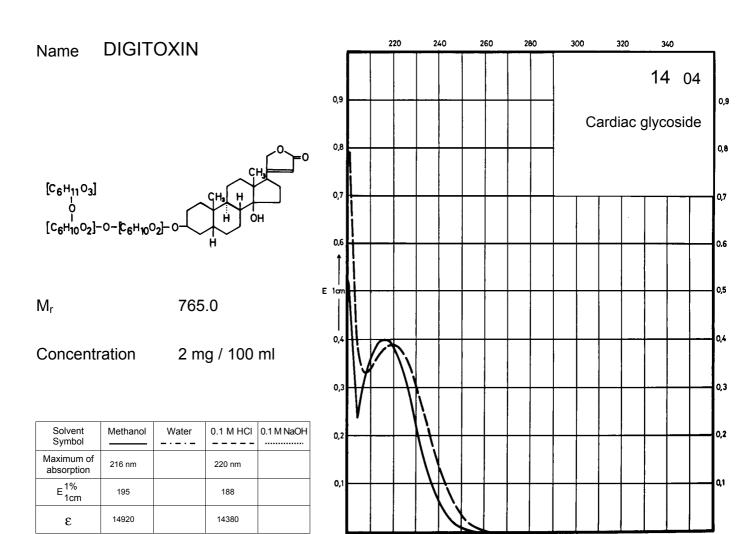

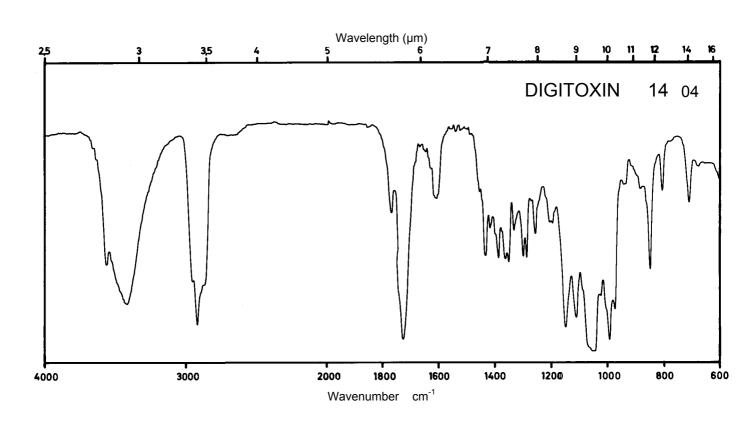





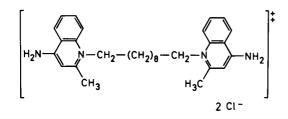

 M_r 527.6


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	338 nm 329 nm 242 nm		327 nm 239 nm	327 nm 239 nm
E ^{1%} 1cm	507 560 862		491 828	444 782
ε	26750 29500 45500		25900 43700	23430 41300

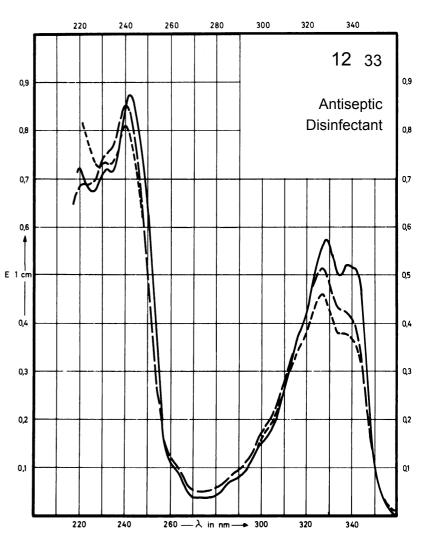

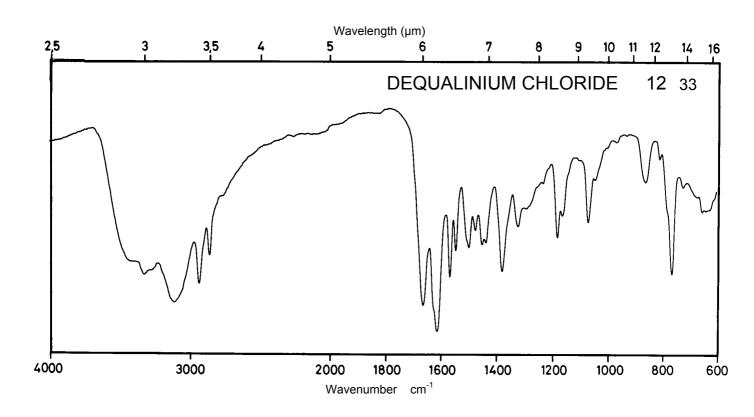

Name 3,5-DIBROMO-4-HYDROXY-

BENZENESULFONIC ACID SODIUM SALT

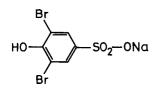



 M_r

354.0


Concentration 3 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 283 nm		289 nm 281 nm	304 nm 258 nm
E ^{1%} _{1cm}	56 54		41 42	134 291
ε	1980 1900		1460 1500	4750 10300

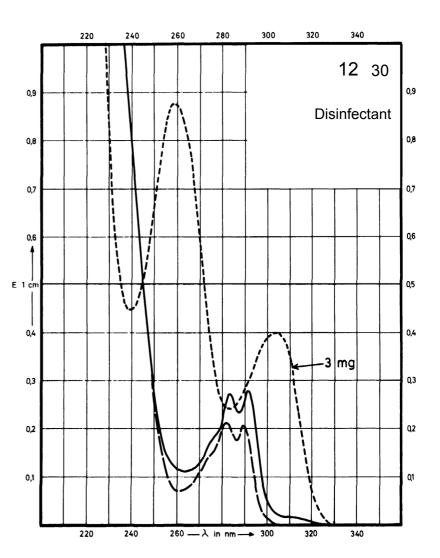

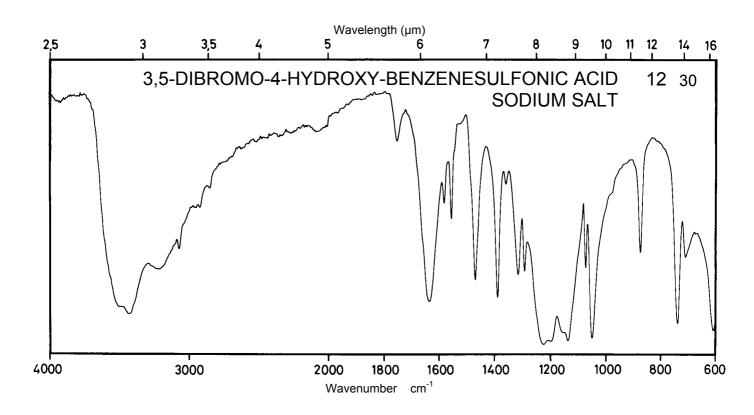



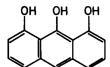



226.2

Concentration

1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	360 nm 288 nm 256 nm			276 nm
E ^{1%} 1cm	418 358 467			490
ε	9460 8100 10560			11080

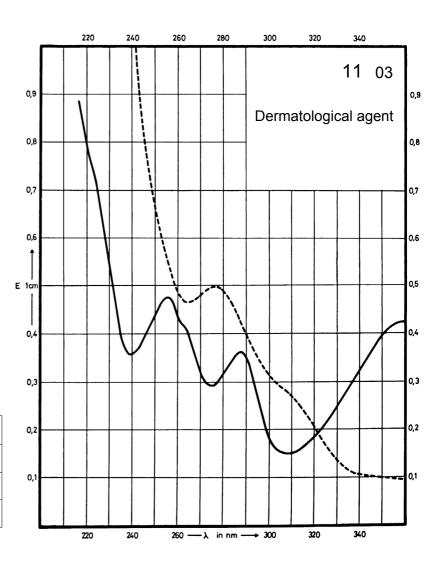

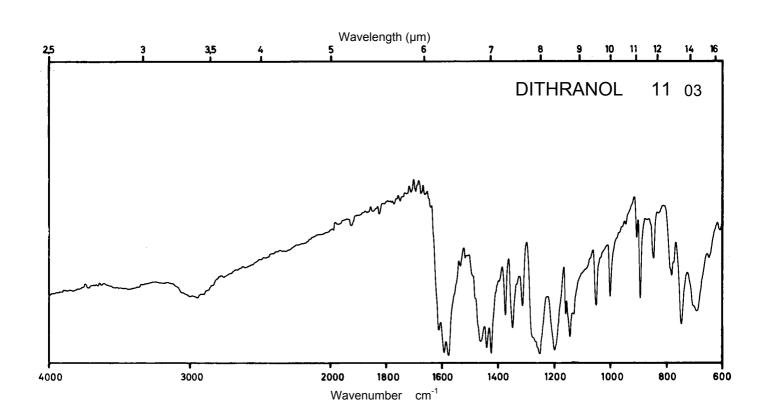



Name DIACETYLAMINO-AZOTOLUENE

 M_{r}


309.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	327 nm			
E ^{1%} 1cm	658			
3	20360			

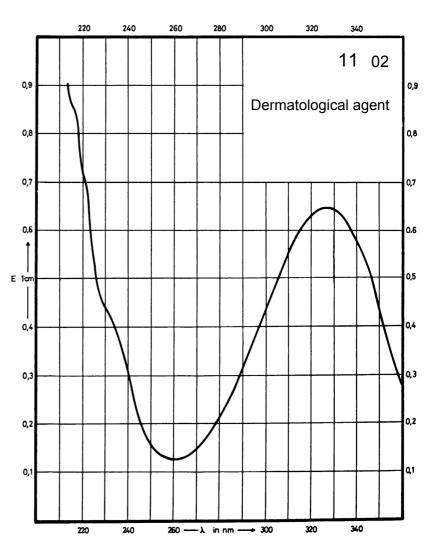

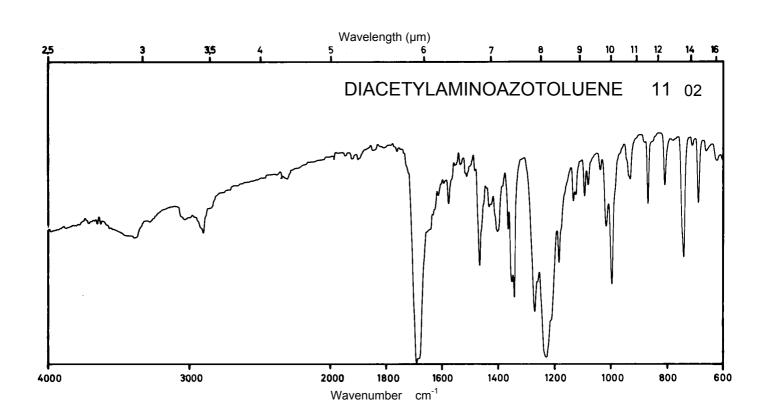





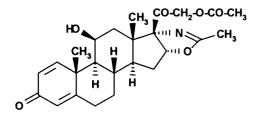

441.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm	246 nm	246 nm	
E ^{1%} 1cm	349	350	355	
ε	15400	15400	15600	

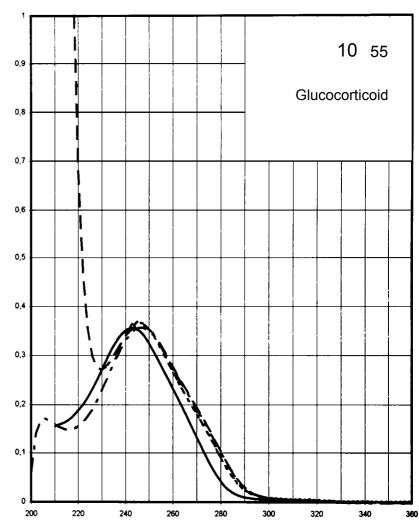

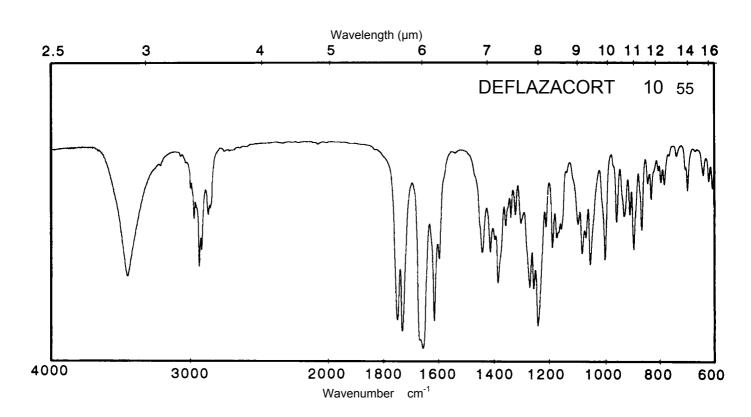



Name DIFLUCORTOLONE 21-VALERATE

 M_{r}


478.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	237 nm			
E ^{1%} 1cm	347			
ε	16620			

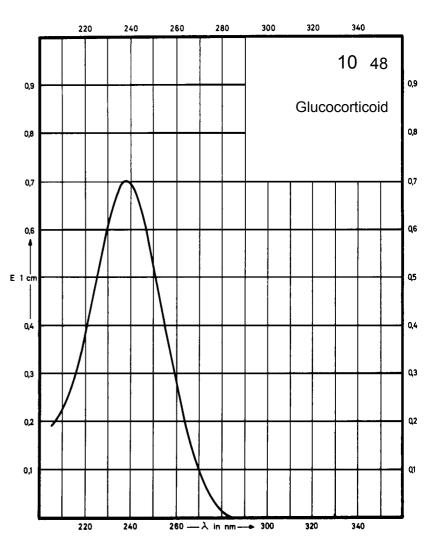

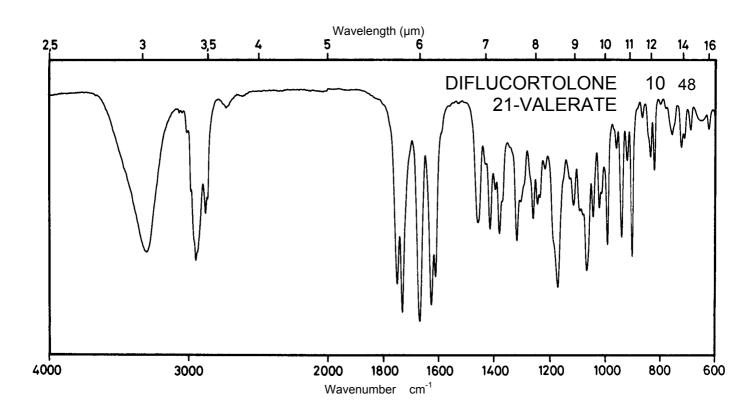





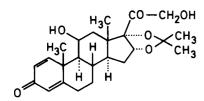


416.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm			
E ^{1%} 1cm	363			
3	15110			

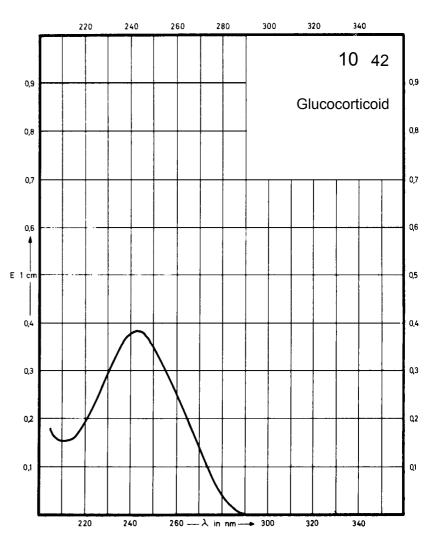

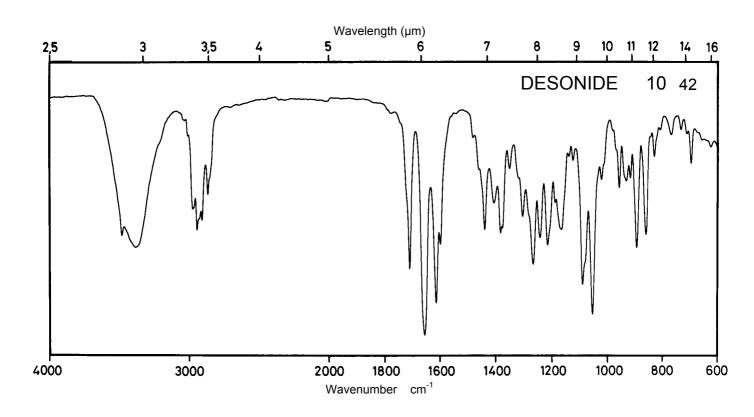

Name DEXAMETHASONE 21-(3,6,9-TRIOXA)-UNDECANOATE

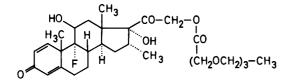


 M_{r}


566.7

Concentration


3.3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm	241 nm		
E ^{1%} _{1cm}	268	271		
ε	15200	15400		

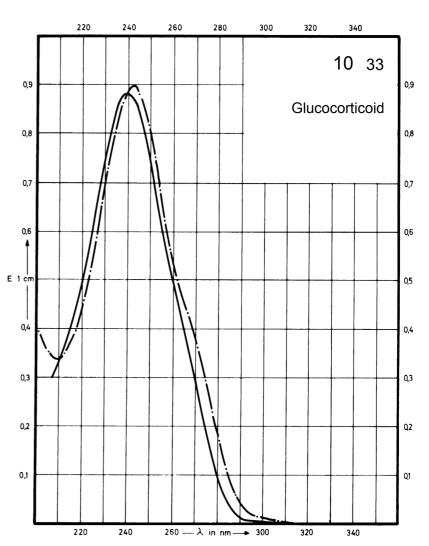

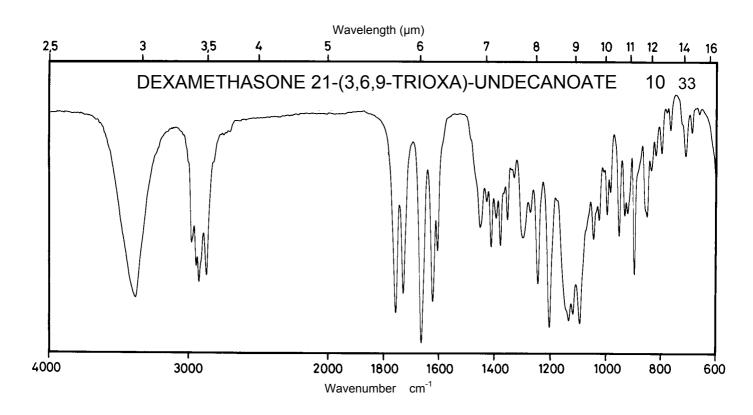


Name DIFLORASONE 17α, 21-DIACETATE

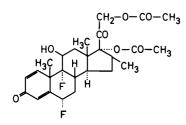


 $\,M_r\,$


494.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm	241 nm		
E ^{1%} 1cm	343	347		
ε	17000	17200		

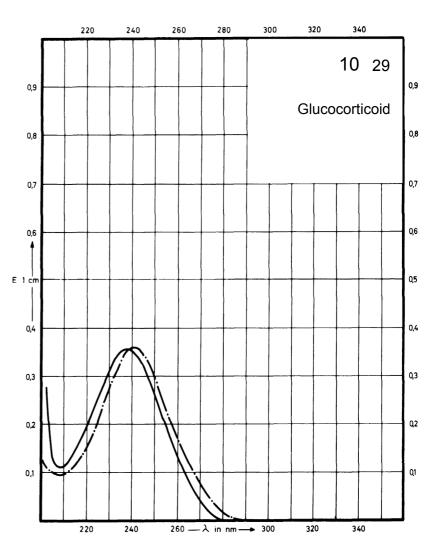

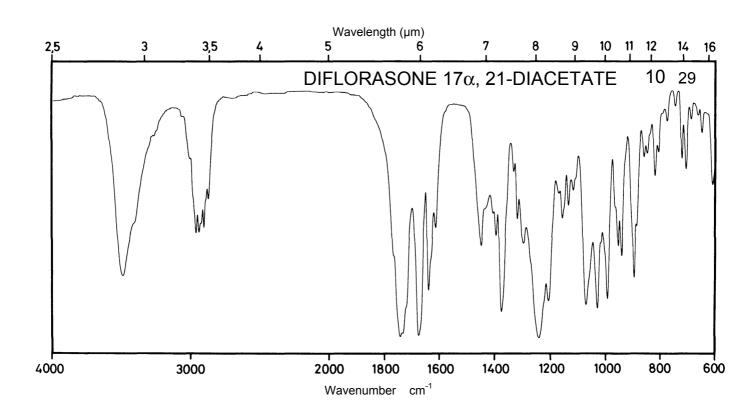





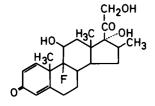


 M_r 392.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm		242 nm	242 nm
E ^{1%} 1cm	386		389	389
ε	15150		15270	15270

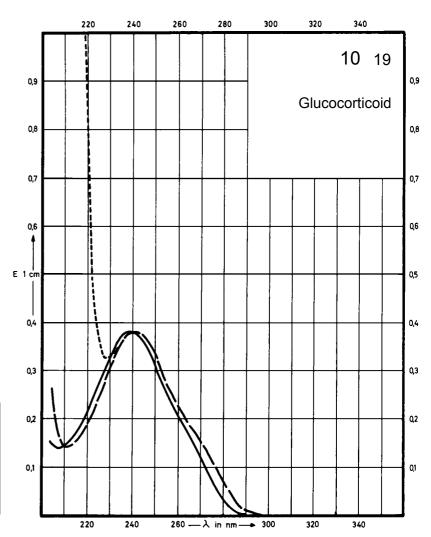

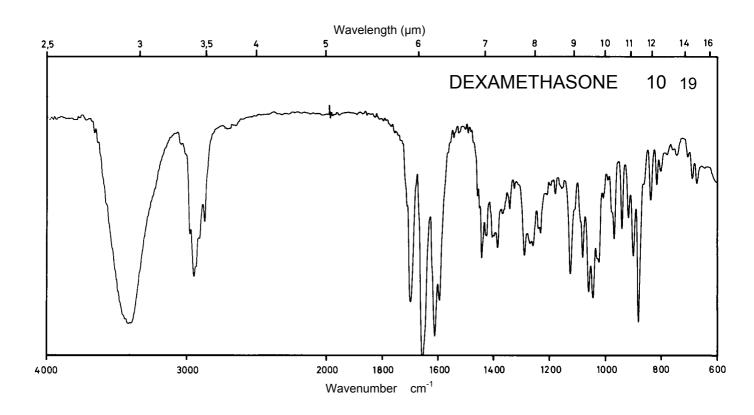





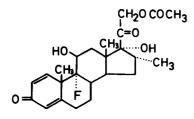


434.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm			
E ^{1%} 1cm	365			
ε	15860			

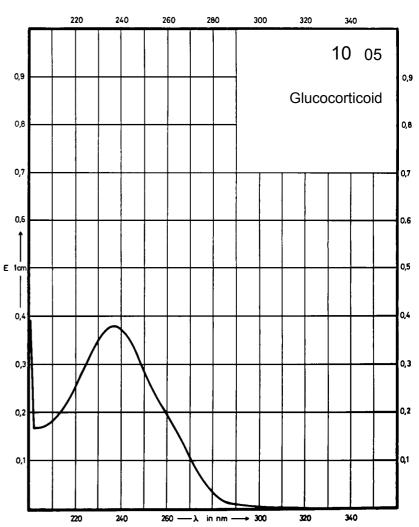

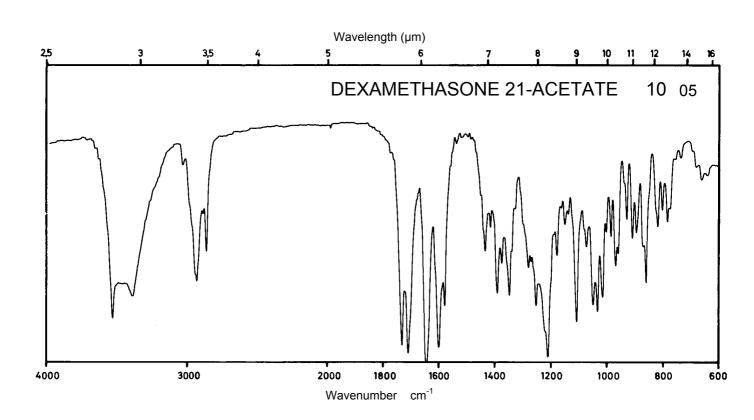





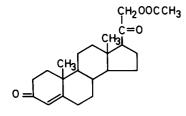


372.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm			
E ^{1%} 1cm	444			
ε	16540			

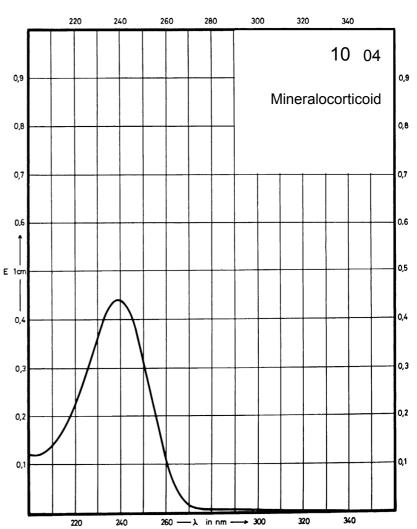

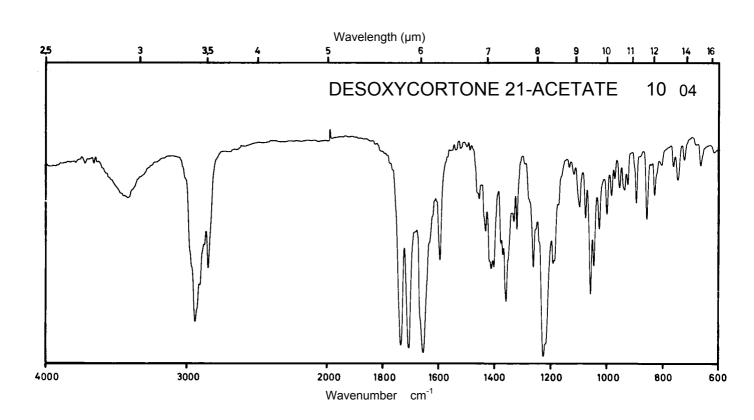





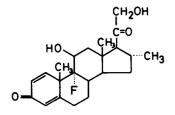


 M_r 376.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm			
E ^{1%} 1cm	429			
ε	16150			

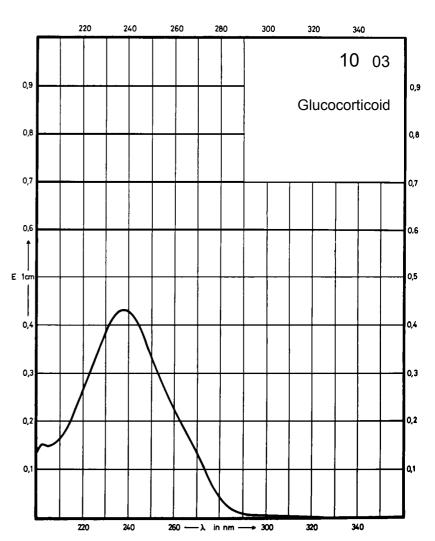

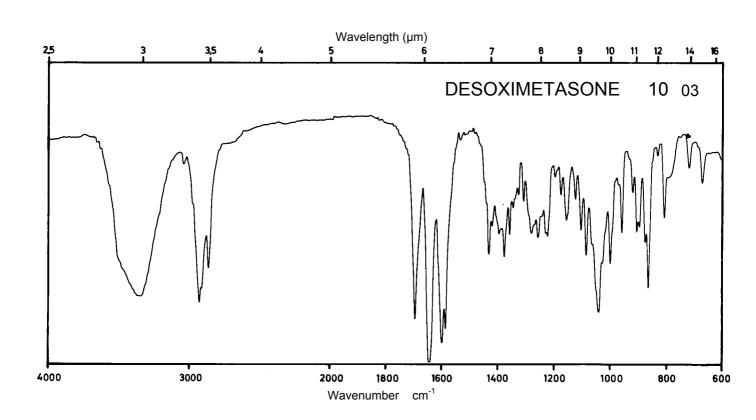



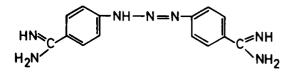



+ 2 CH3CONH CH2 COOH + 4 H2O

 M_{r}


587.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	375 nm		257 nm 224 nm	425 nm 247 nm
E ^{1%} 1cm	635		290 278	666 239
3	37310		17040 16340	39130 14040

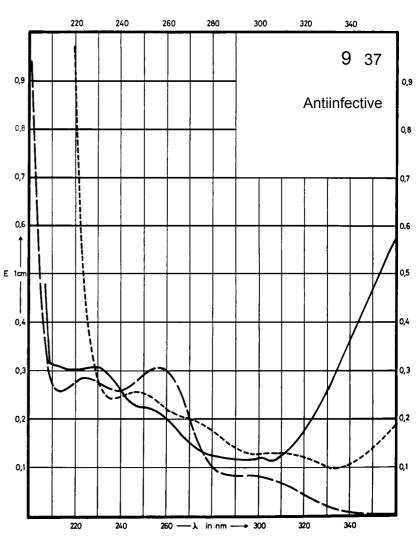

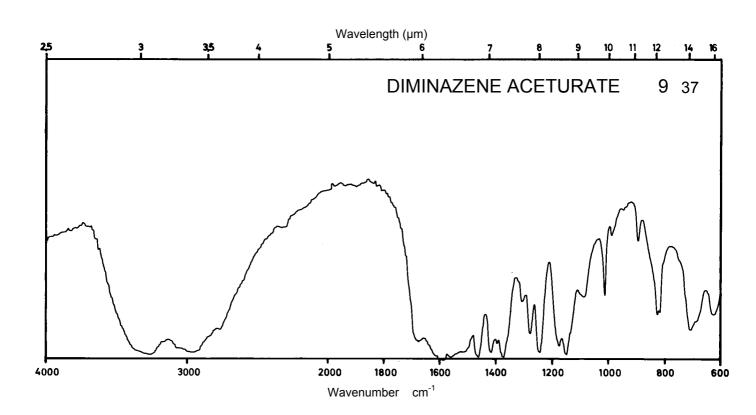





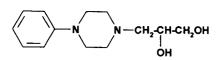



236.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm	234 nm		
E ^{1%} _{1cm}	508	398		
ε	12000	9400		

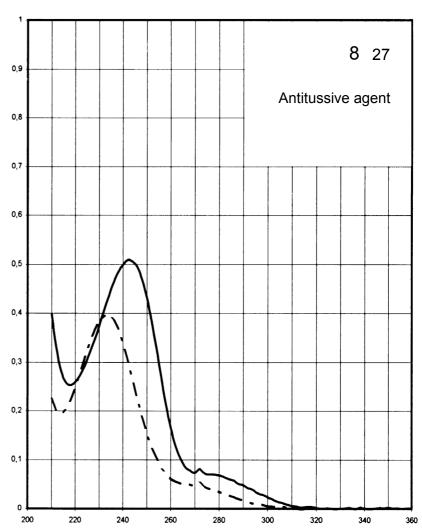

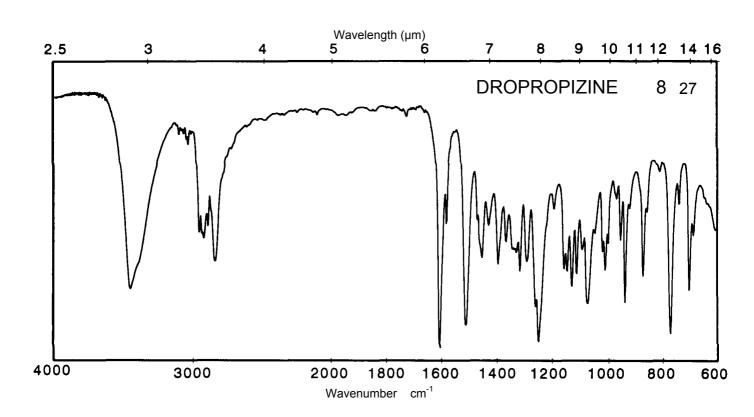

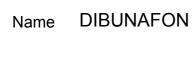

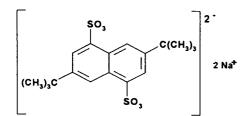


 M_r


444.5

Concentration (


0.25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	326 nm	328 nm	328 nm	328 nm
	283 nm	282 nm	283 nm	282 nm
	234 nm	234 nm	234 nm	234 nm
E ^{1%} 1cm	74	91	93	93
	161	155	158	156
	2070	1990	1980	1960
ε	3270	4040	4150	4100
	7140	6880	7000	6900
	91900	88400	88100	87300

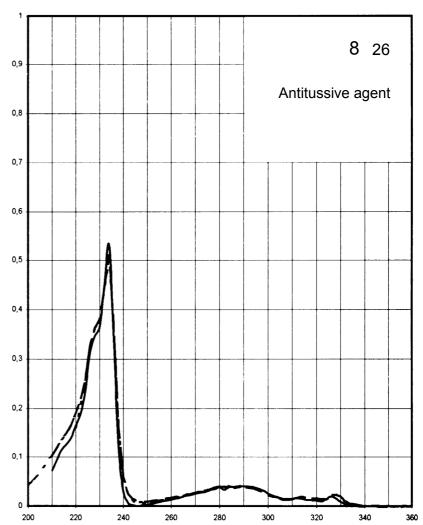

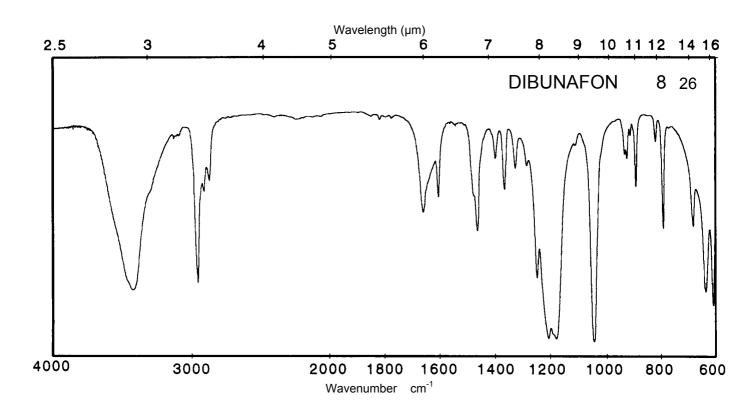


Name DEXTROMETORPHAN HYDROBROMIDE

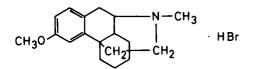


 M_{r}


352.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm 280 nm		278 nm	279 nm
E ^{1%} 1cm	55 60		53	57
ε	1940 2110		1870	2010

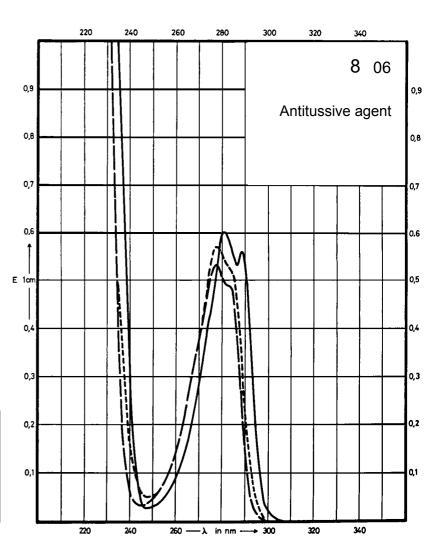

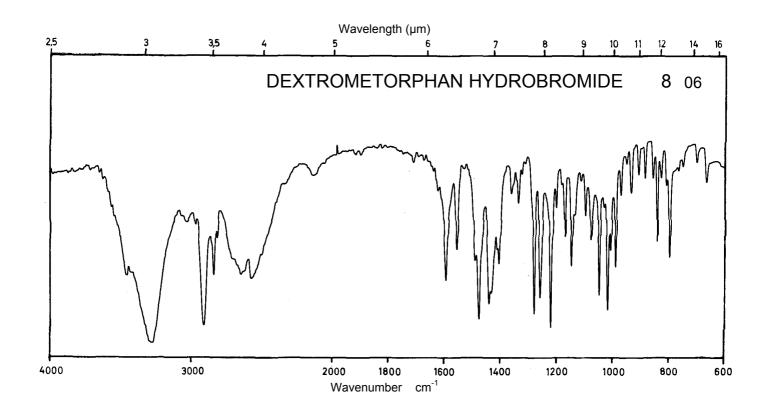





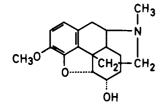

 M_r 301.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm			
E ^{1%} 1cm	53			
ε	1600			

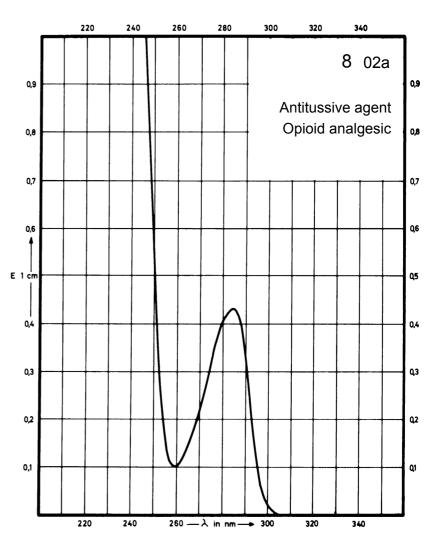

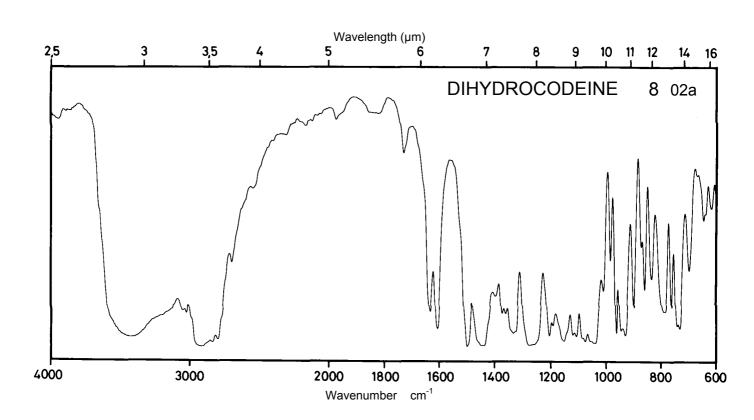


Name DIHYDROCODEINE TARTRATE


$$H_3$$
C-O CH_2 CH_2 CH_2 $COOH$ $COOH$ $COOH$

 M_{r}


451.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm		283 nm	282 nm
E ^{1%} 1cm	39.8		37.6	35.0
ε	1800		1700	1580

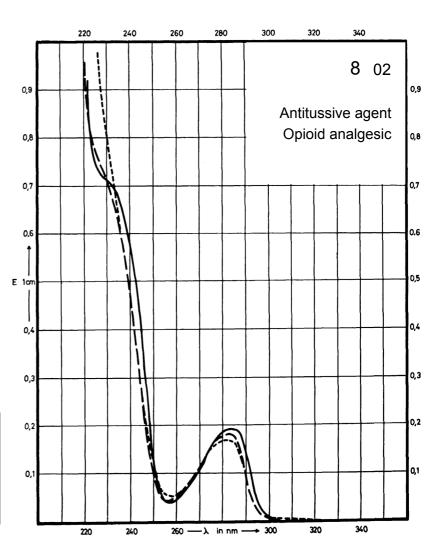

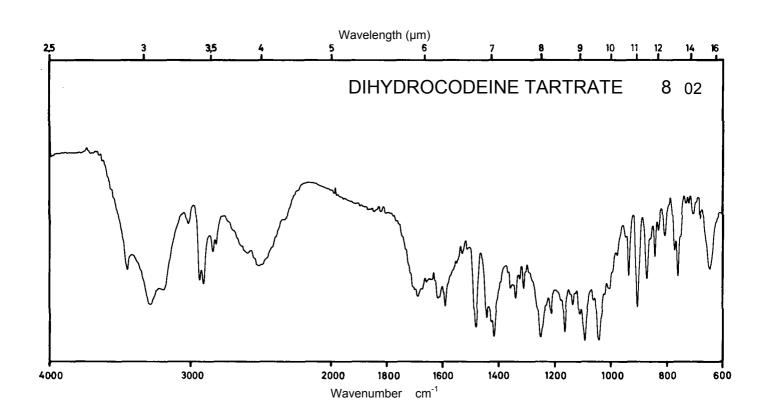

Name DICLOFENAC SODIUM SALT

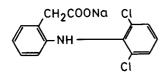


 M_{r}


282.7

Concentration 1 m


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm		274 nm	275 nm
E ^{1%} _{1cm}	425		288	327
3	12030		8150	9250

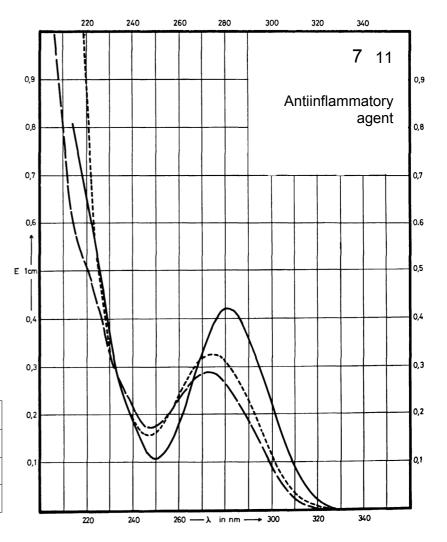

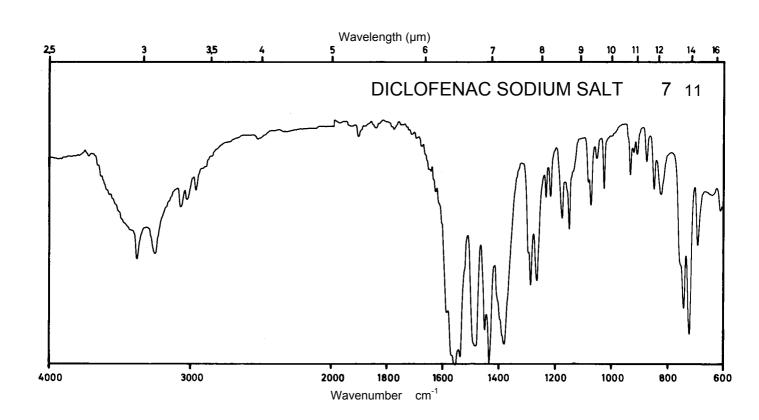



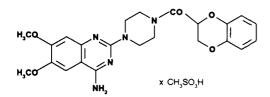



547.6

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	341 nm 330 nm 246 nm	330 nm 246 nm	330 nm 246 nm	345 nm 275 nm 251 nm
E ^{1%} 1cm	183 197 1015	179 1021	199 1040	100 343 997
ε	1000 10800 55600	9800 55900	10900 57000	5500 18800 54600

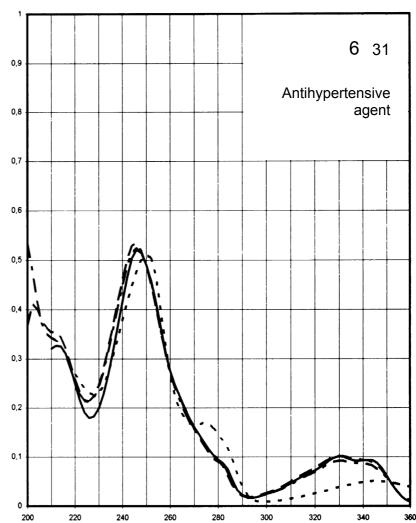

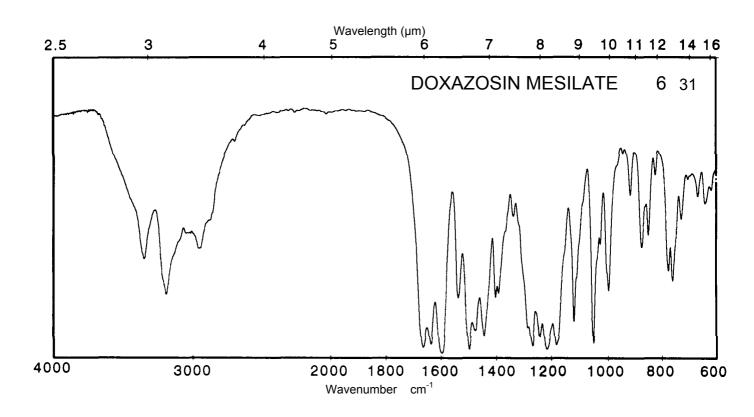




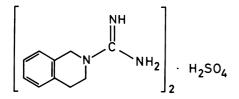


 M_r 448.5


Concentration 2 mg / 100 ml 40 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm	262 nm	262 nm	
E ^{1%} 1cm	13.2	13.7	13.8	
ε	593	614	617	

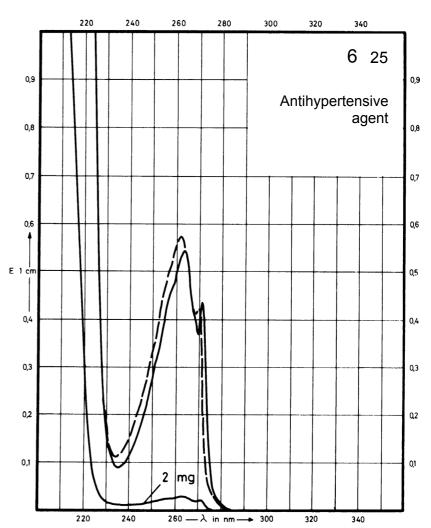

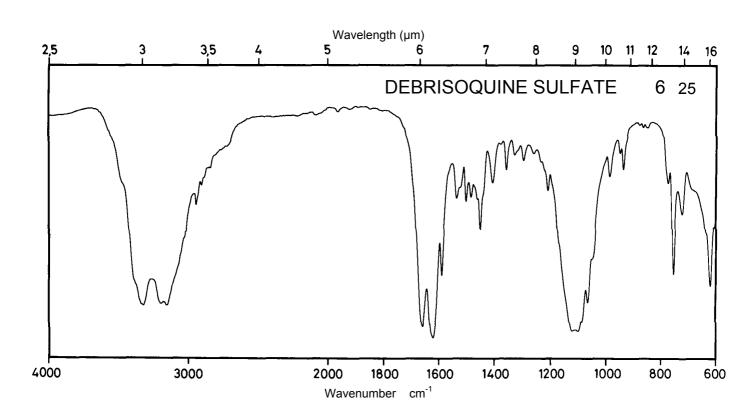





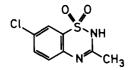


 M_r 230.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			264 nm	280 nm
E ^{1%} 1cm			502	580
ε			11600	13400

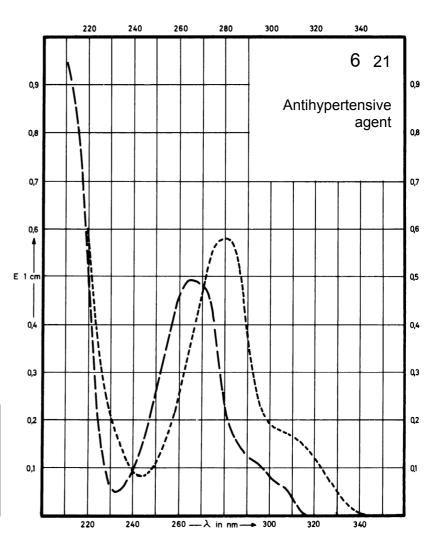

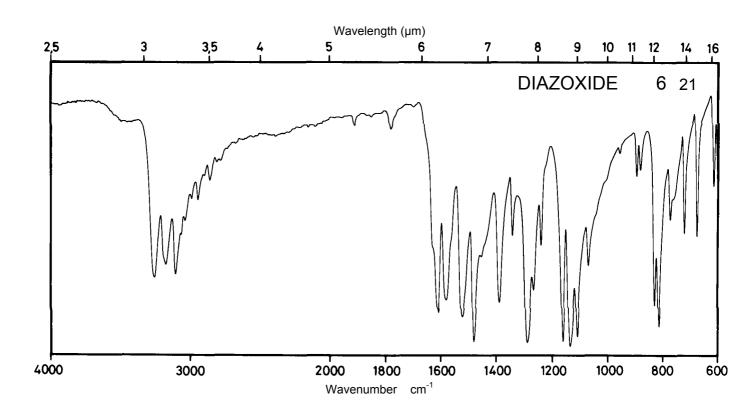





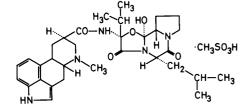


 M_r 673.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 280 nm		279 nm 220 nm	279 nm
E ^{1%} 1cm	86 102		95 482	93
ε	5800 6800		6400 32500	6200

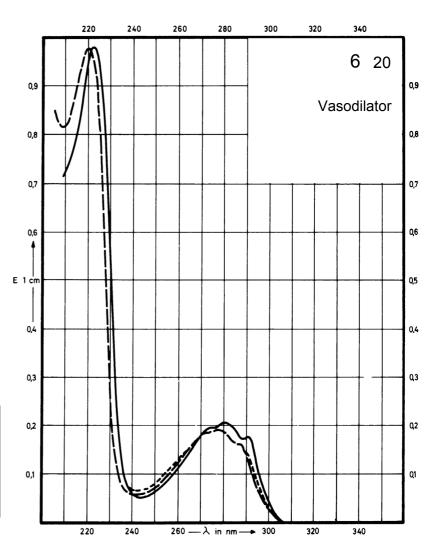

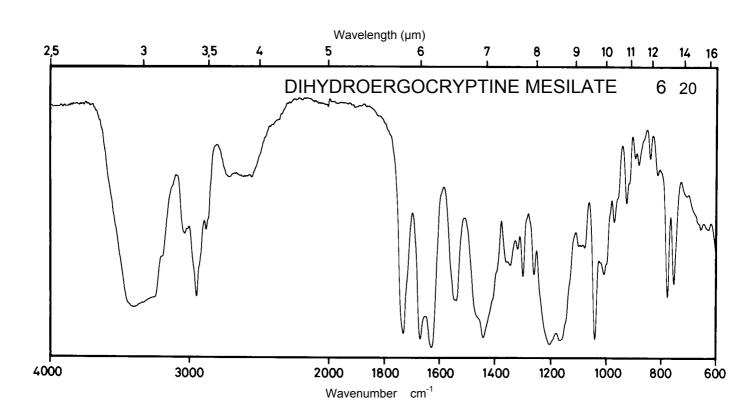





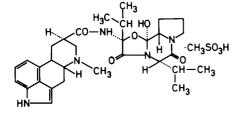


659.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 280 nm		279 nm 220 nm	279 nm
E ^{1%} 1cm	87 103		96 495	93
ε	5800 6800		6400 32700	6200

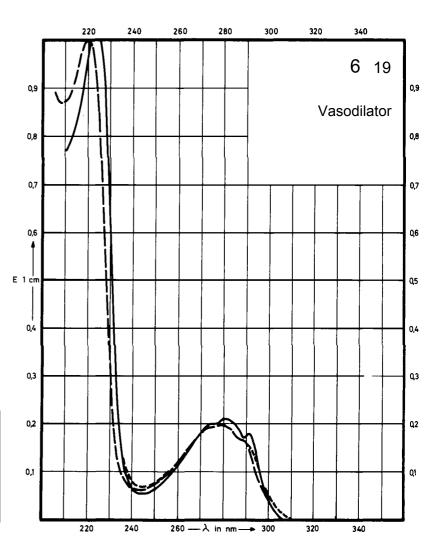

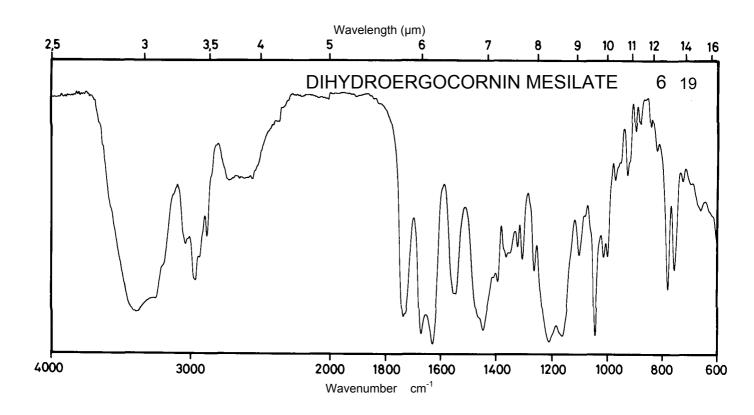

Name DIHYDRALAZINE HYDROGEN SULFATE

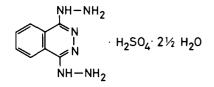



 M_r

333.3


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			306 nm 274 nm 240 nm	
E ^{1%} *			213 204 377	
ε			6130 5870 10870	

^{*} Calculated on dried substance

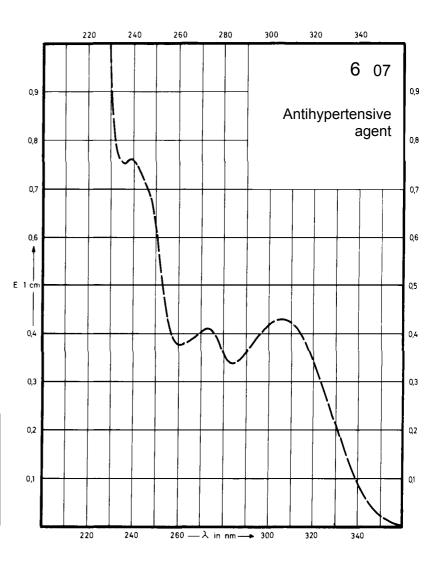

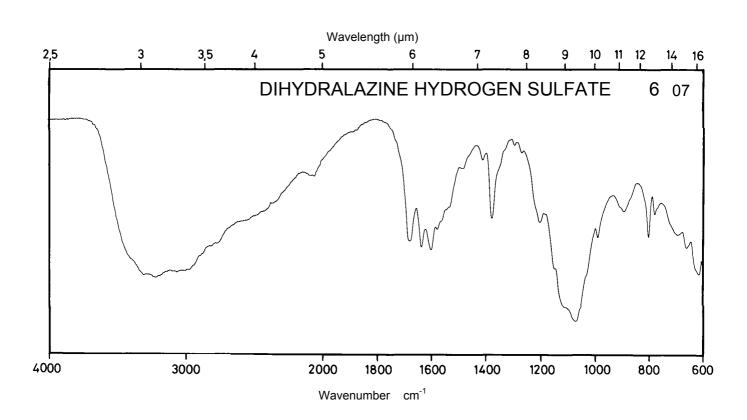





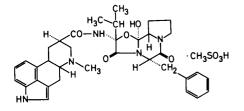

 M_r 707.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 280 nm		279 nm 220 nm	279 nm
E ^{1%} 1cm	83 94		90 517	88
ε	5900 7000		6400 36600	6200

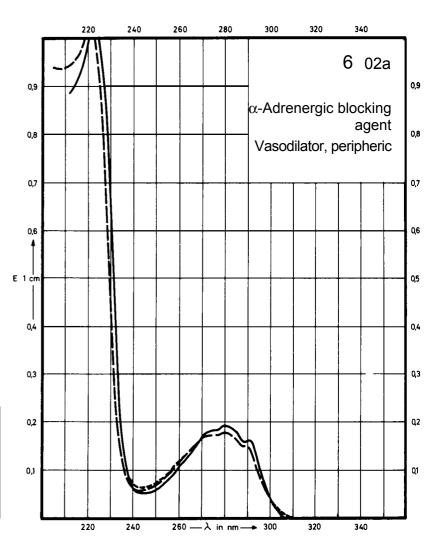

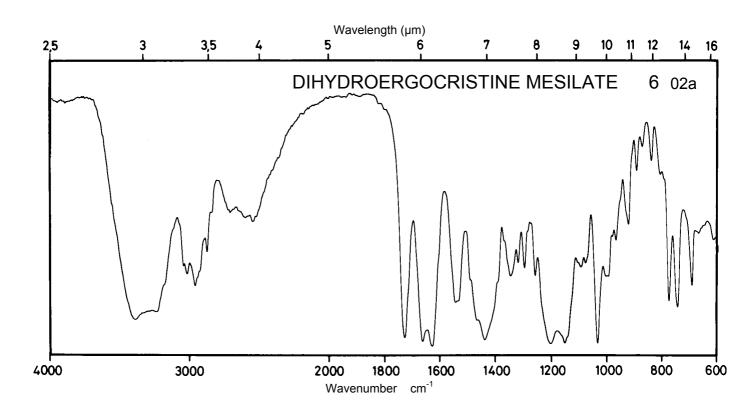


Name DIHYDROERGO-CRISTINE


$$\begin{array}{c} & & & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ &$$

 M_r


611.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		279 nm	280 nm
E ^{1%} 1cm	111		103	100
ε	6790		6300	6120

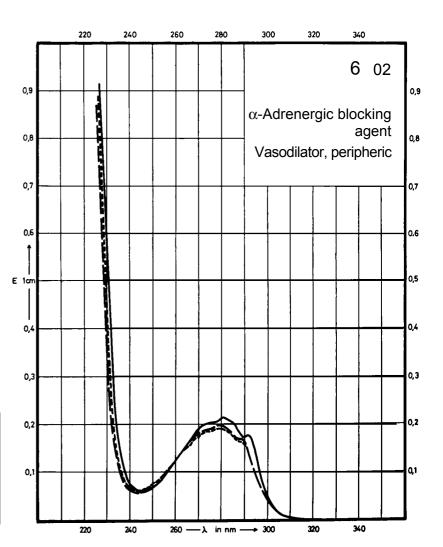

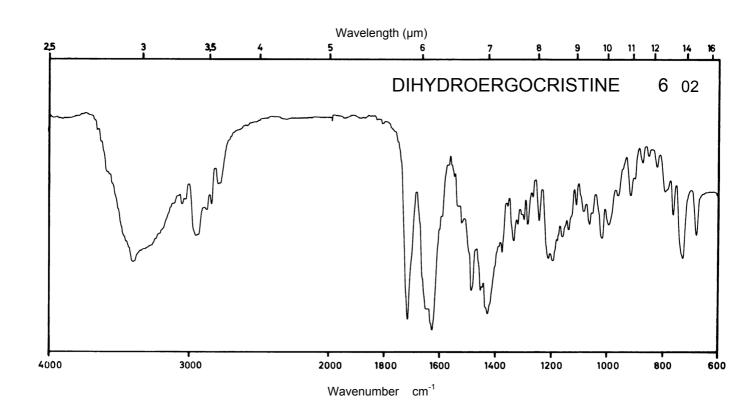


Name DOXYCYCLINE HYDROCHLORIDE

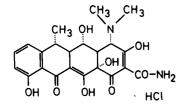


 M_{r}

480.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	352 nm 269 nm		346 nm 269 nm	376 nm 255 nm
E ^{1%} 1cm	297 368		246 381	329 326
ε	14290 17700		11850 18330	15820 15680

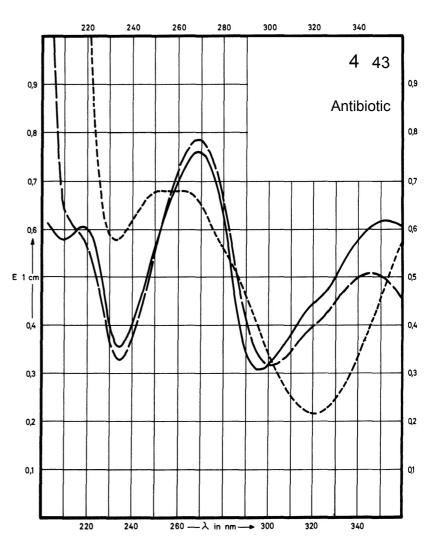

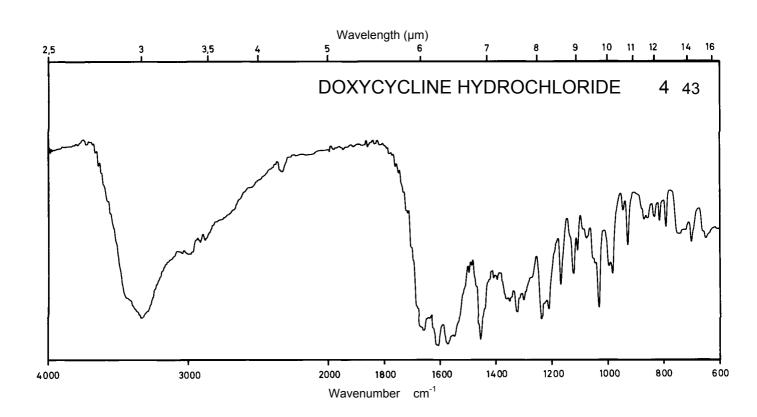



Name DICLOXACILLIN SODIUM SALT

 M_{r} 510.3


Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 275 nm	282 nm 275 nm		281 nm 274 nm
E ^{1%} 1cm	11.7 12.3	13.3 13.3		13.3 14.0
ε	596 630	676 676		676 715

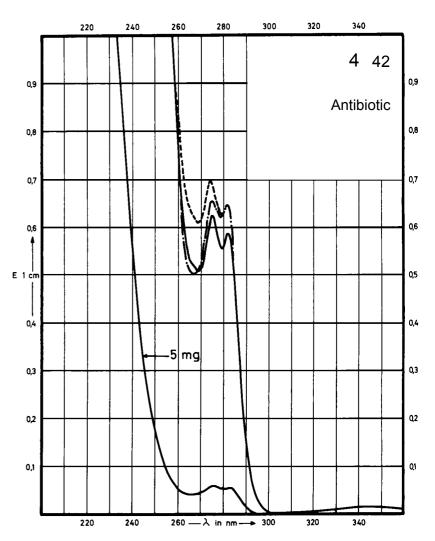

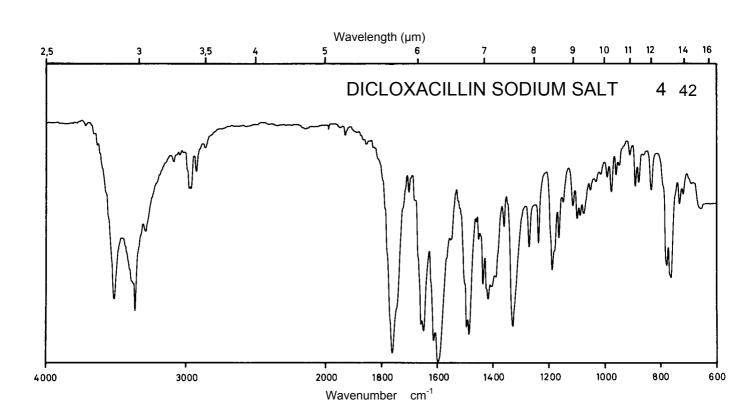


Name DEMECLOCYCLINE HYDROCHLORIDE

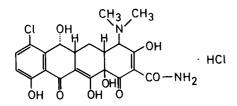



 M_r


501.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	373 nm 268 nm 229 nm		367 nm 267 nm 227 nm	388 nm 276 nm 241 nm
E ^{1%} 1cm	306 317 384		252 364 356	347 282 338
3	15330 15890 19270		12600 18250 17840	17380 14150 16970

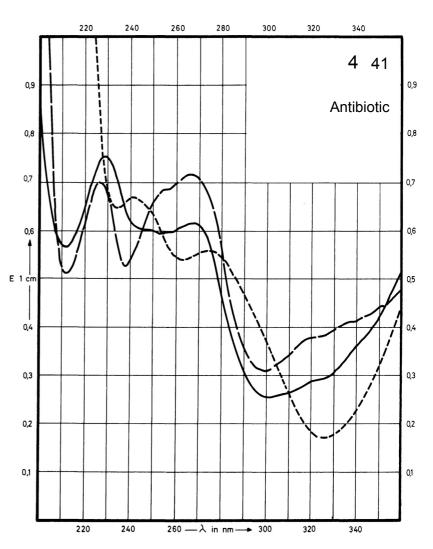


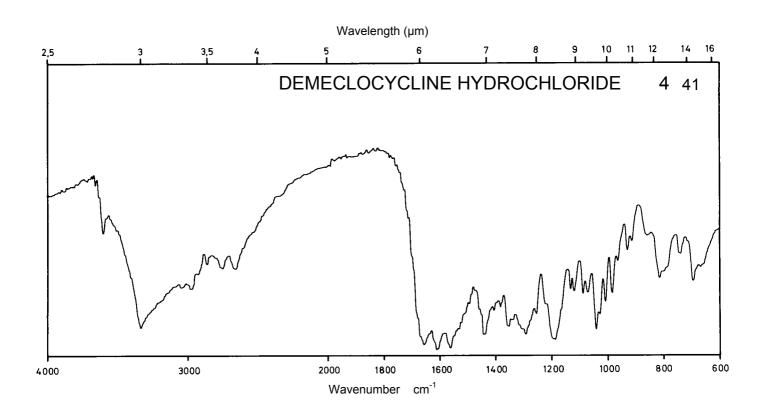

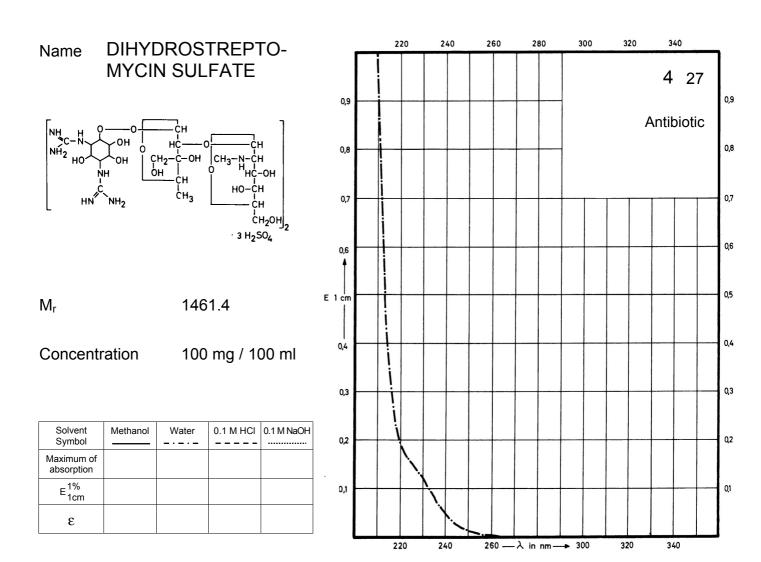

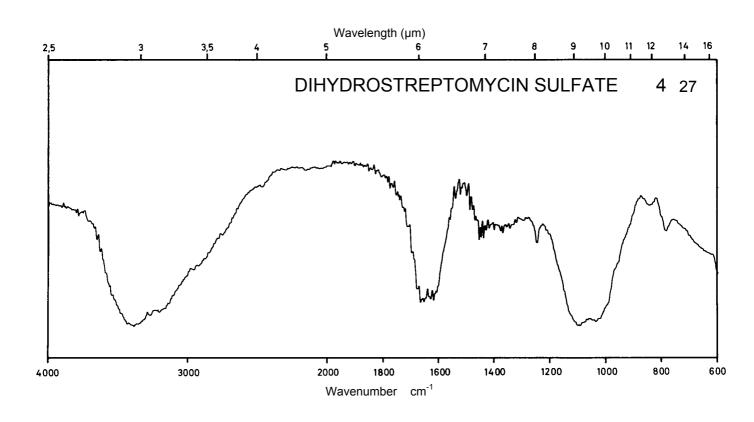





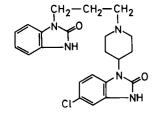


 M_r 425.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm 230 nm		284 nm 228 nm	292 nm 245 nm
E ^{1%} 1cm	294 270		286 295	322 222
ε	12500 11500		12200 12500	13700 9400

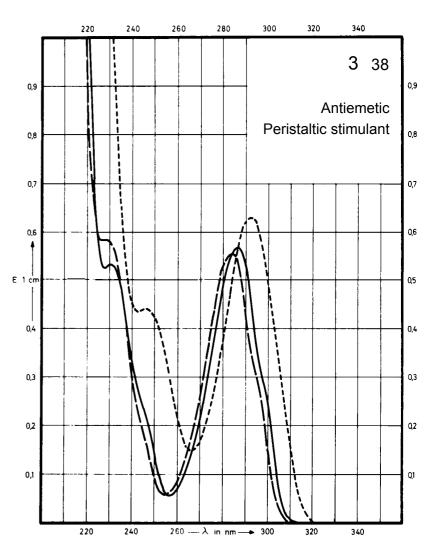

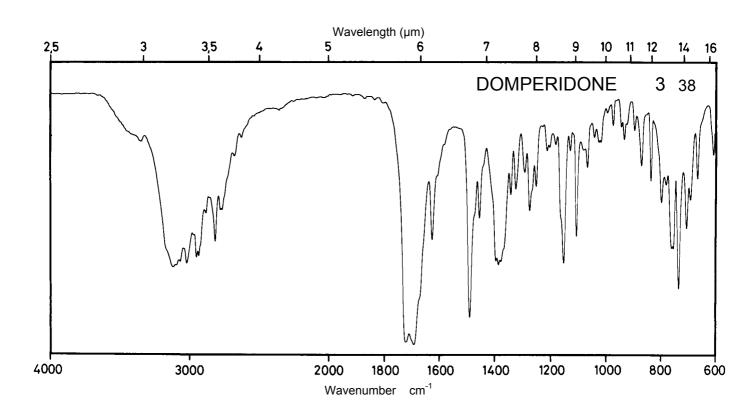




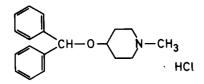


 M_r


317.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm		258 nm	
E ^{1%} 1cm	13.0		14.3	
3	410		450	

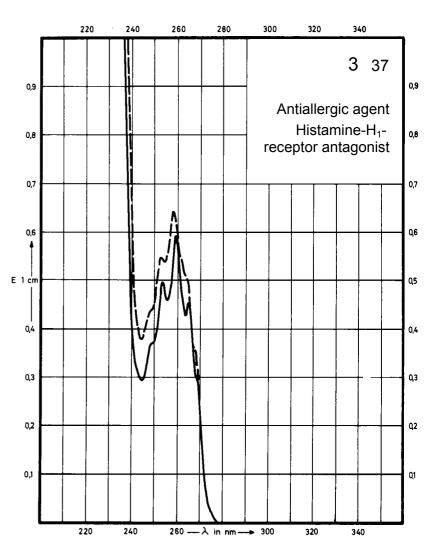

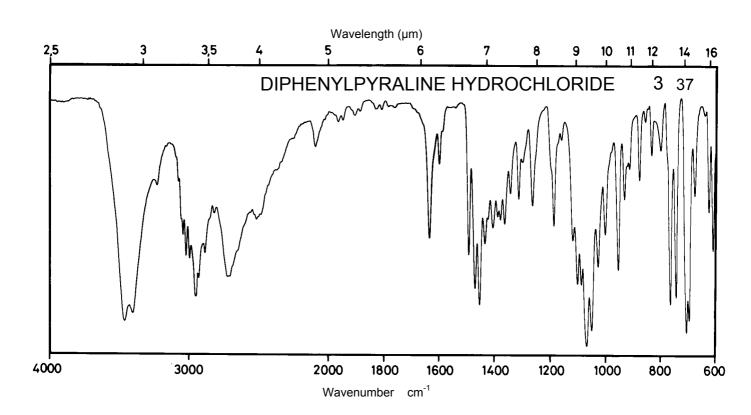



Name DIFENIDOL HYDROCHLORIDE

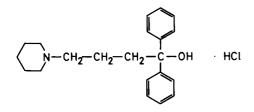


 M_r 345.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm		257 nm	
E ^{1%} 1cm	12.8		13.3	
ε	440		460	

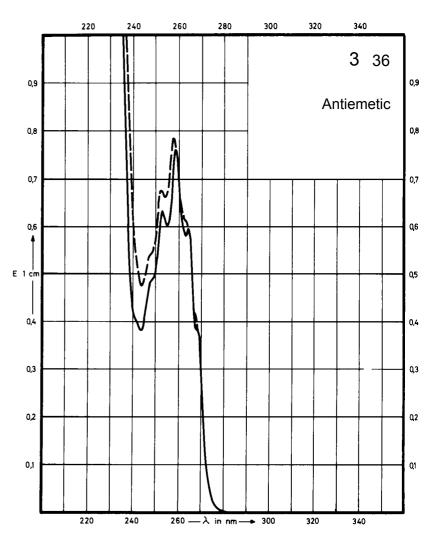

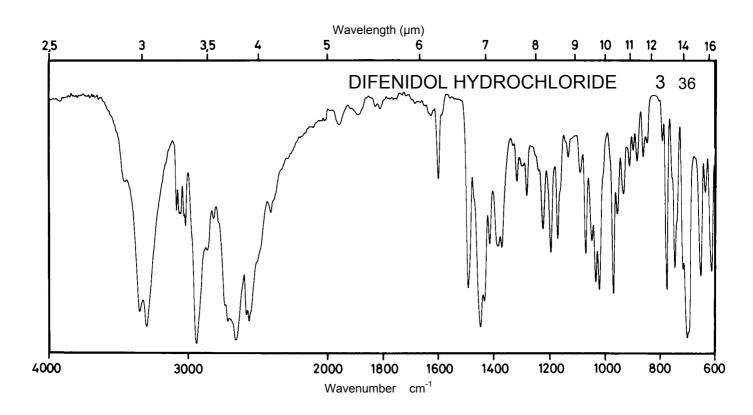





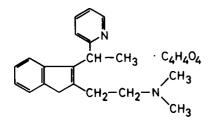


 M_{r} 408.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm	258 nm	260 nm	261 nm
E ^{1%} 1cm	370	392	438	404
ε	15120	16030	17900	16520

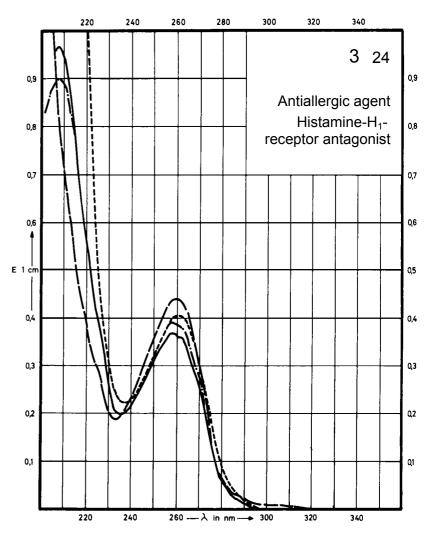

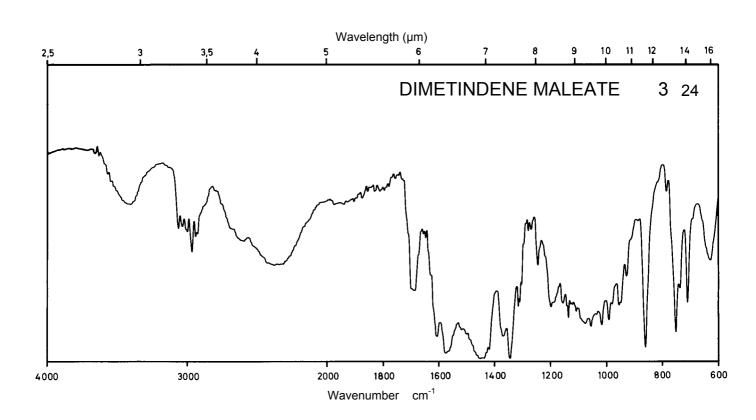





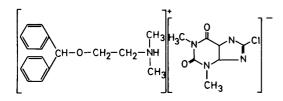


470.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm		277 nm	279 nm
E ^{1%} 1cm	265		261	285
ε	12460		12270	13400

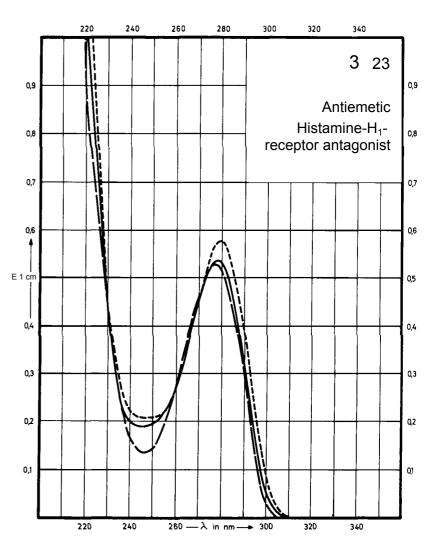

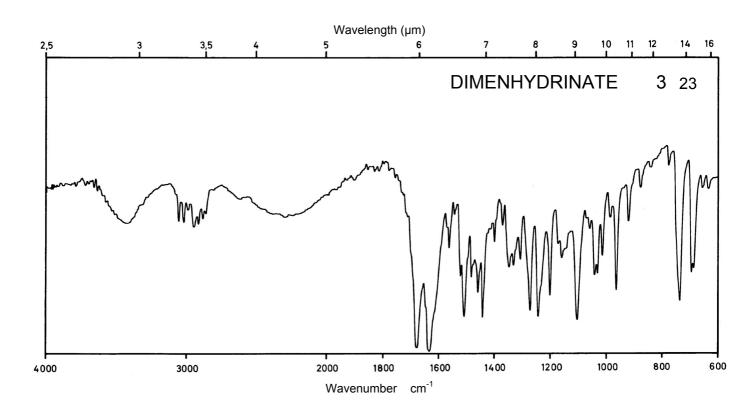


Name TOLPROPAMINE HYDROCHLORIDE

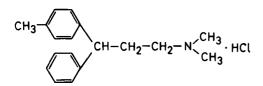


 M_{r}

289.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm 265 nm 274 nm	260 nm 264 nm 273 nm	260 nm 264 nm 273 nm	
E ^{1%} 1cm	19.5 20.9 15.0	21.5 21.6 14.4	21.0 21.3 14.2	
ε	565 606 435	623 626 417	609 617 412	

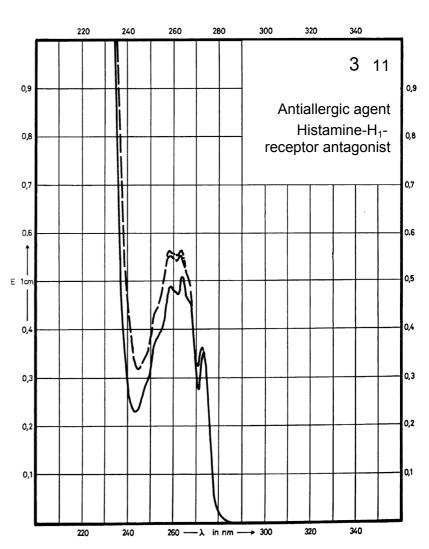

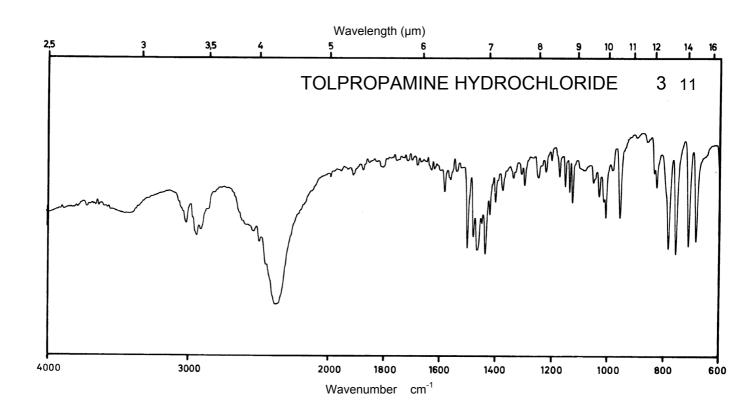



Name DIPHENHYDRAMINE HYDROCHLORIDE

 M_{r}


291.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm		258 nm 252 nm	
E ^{1%} 1cm	12.1 15.5 12.8		13.3 11.4	
ε	353 452 374		388 333	

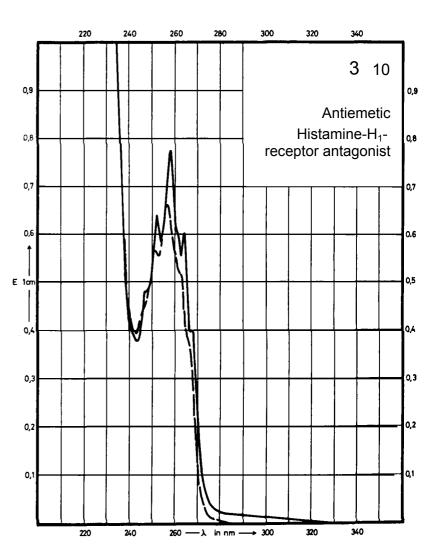

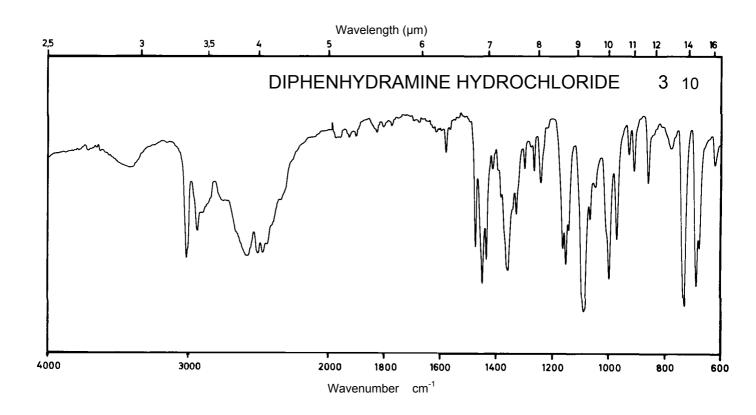





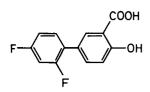


250.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	304 nm 255 nm			273 nm
E ^{1%} 1cm	131 591			453
ε	3280 14880			11300

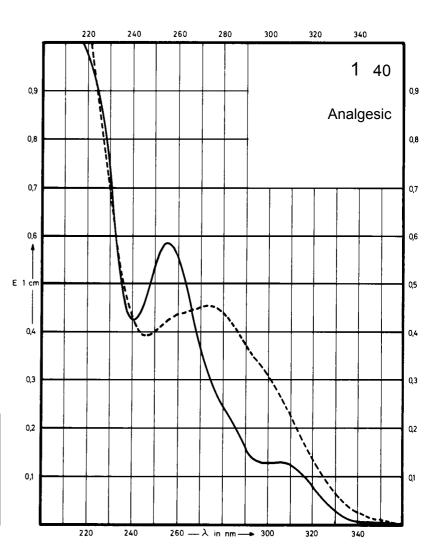

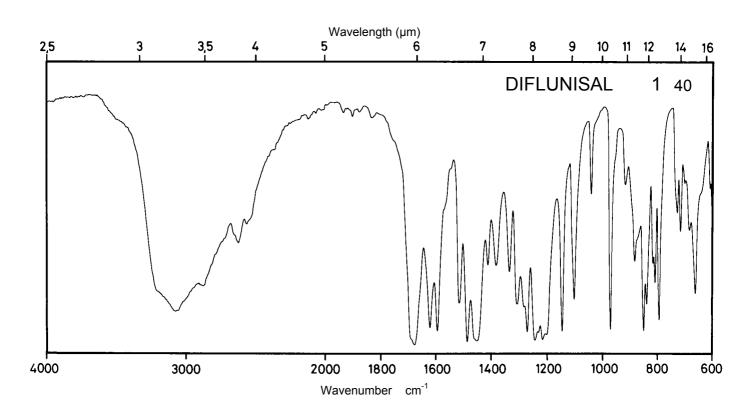


Name DIAMORPHINE HYDROCHLORIDE

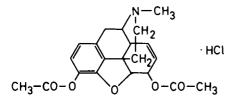


 M_{r}


405.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		279 nm	298 nm
E ^{1%} 1cm	48		42	60
ε	1930		1710	2420

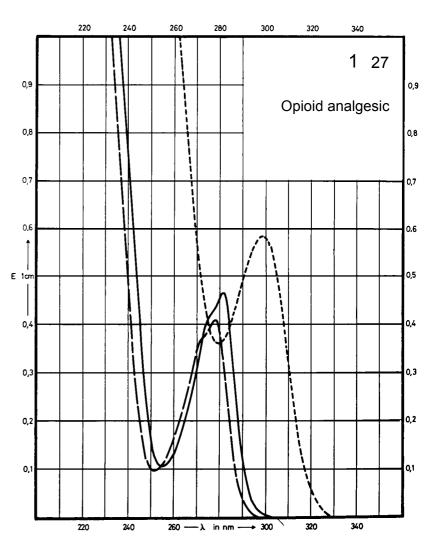

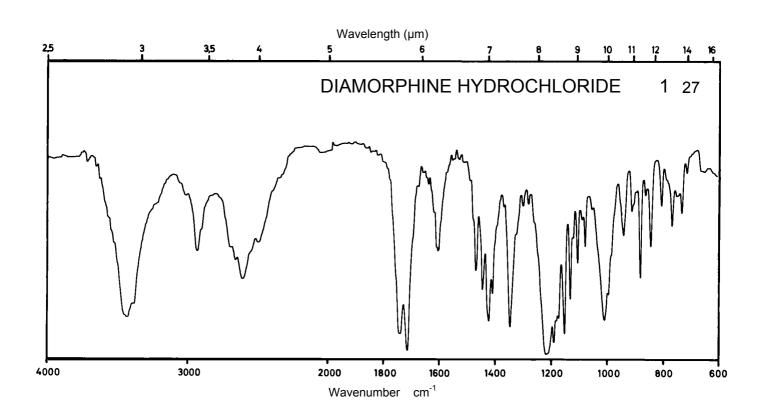





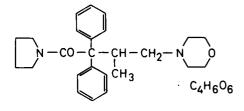


542.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm		260 nm	
E ^{1%} 1cm	7.0		7.5	
ε	380		406	

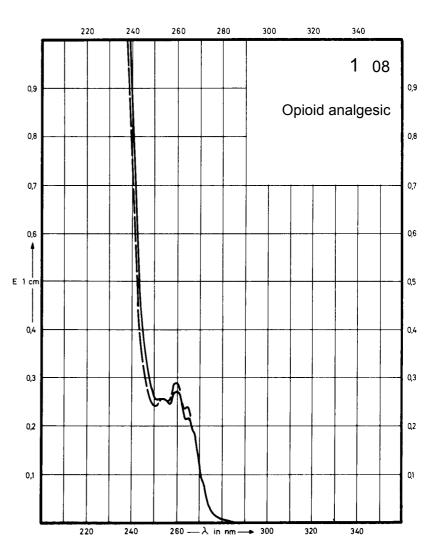

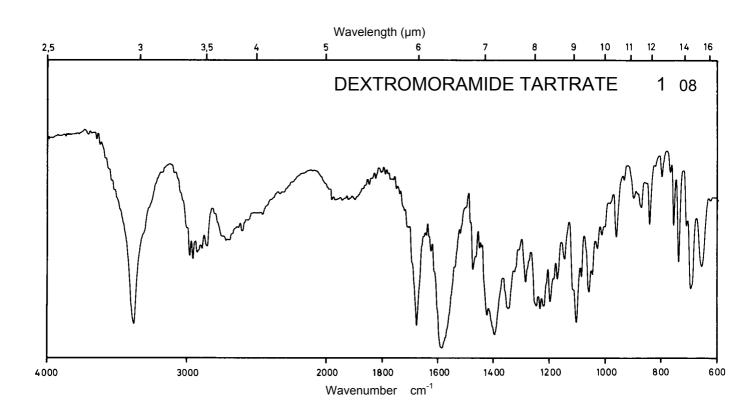


Name DIMETHYLAMINO PHENAZONE

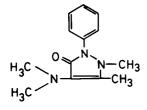


 M_{r}


231.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm 234 nm		256 nm	264 nm
E ^{1%} 1cm	356 408		386	353
ε	8230 9440		8930	8160

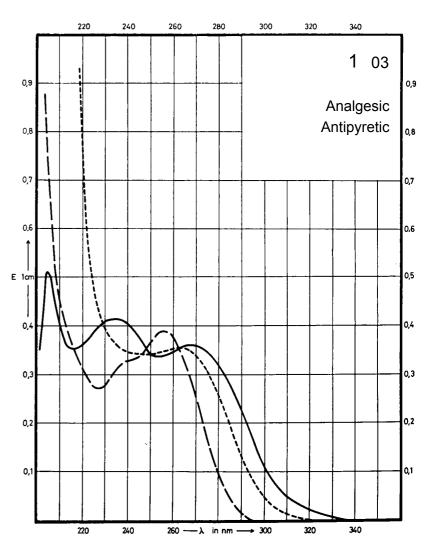

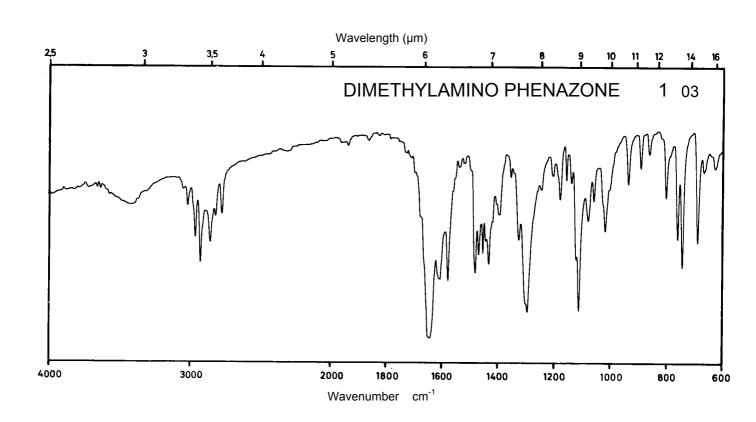




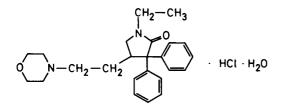



 M_r 433.0


Concentration 2.5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 253 nm	264 nm 258 nm 252 nm		
E ^{1%} 1cm	8.48 9.91 8.48	8.19 10.2 8.59		
ε	367 429 367	355 440 372		

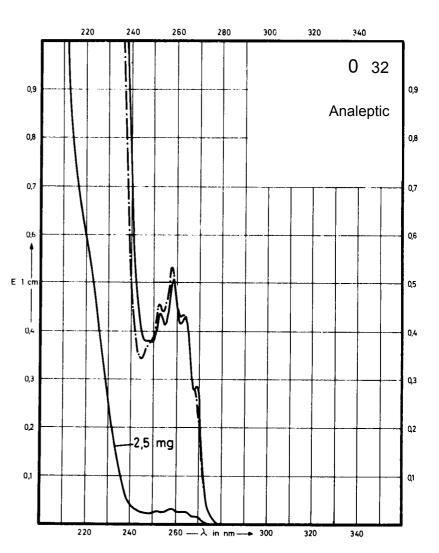


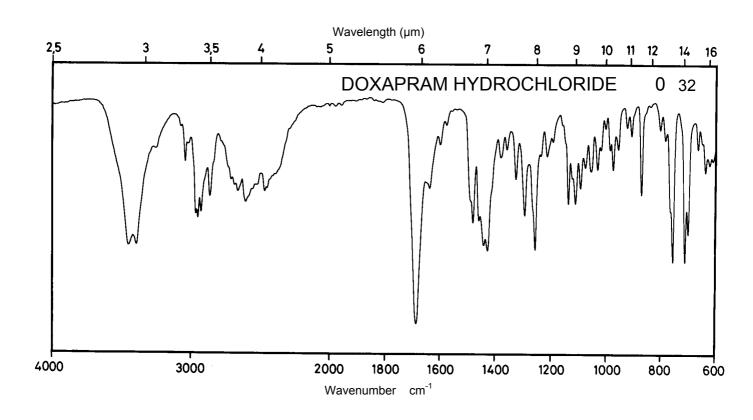

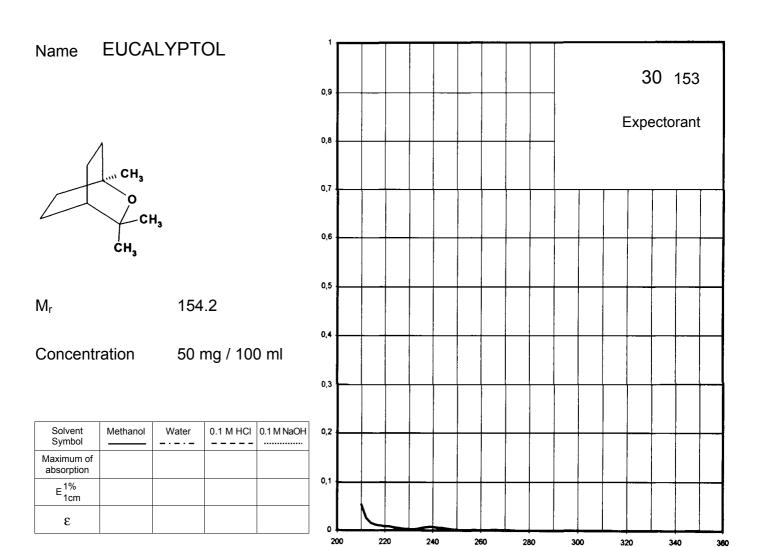

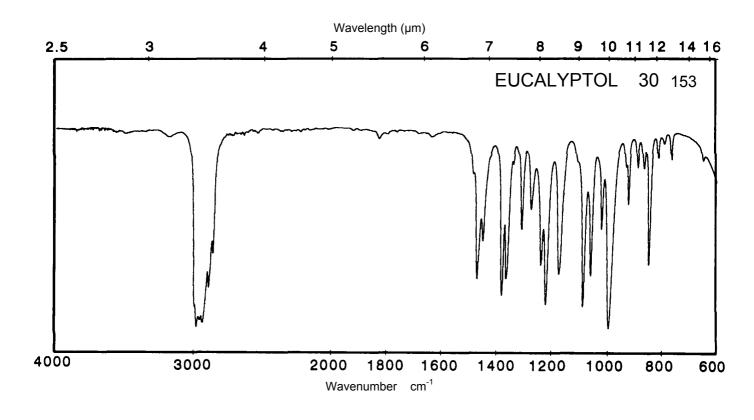




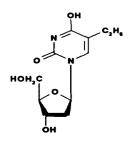


 M_{r}

256.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm	267 nm	267 nm	267 nm
E ^{1%} 1cm	374	380	383	290
3	9600	9750	9800	7400

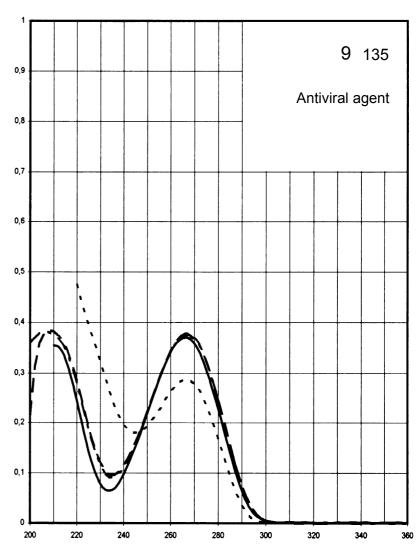

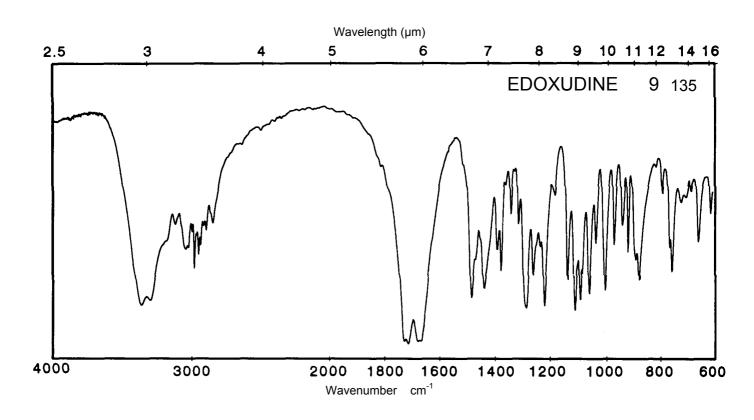


Name ESTRAMUSTINE 17β-(DISODIUM PHOSPHATE)


 M_r

563.8


Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 268 nm			
E ^{1%} 1cm	13.7 14.0			
3	770 790			

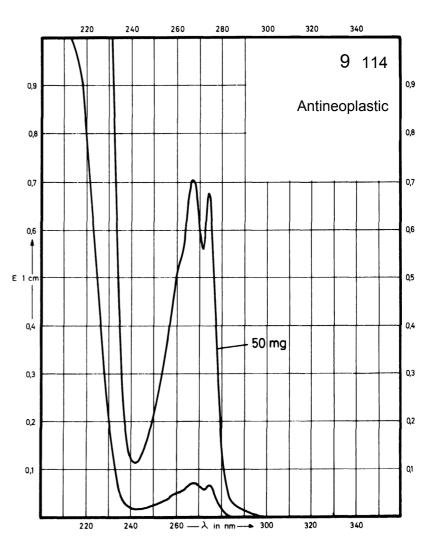

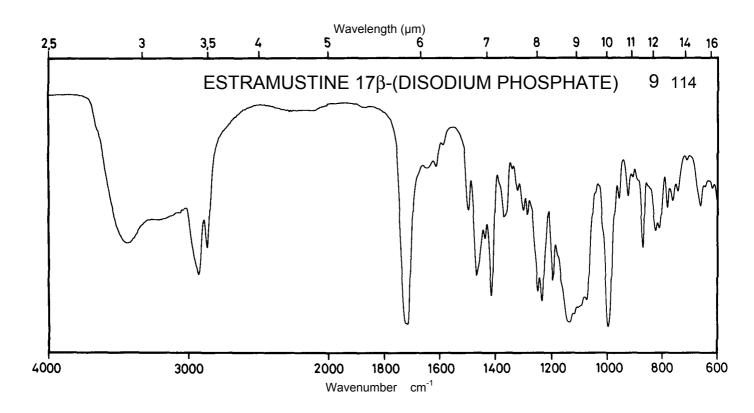





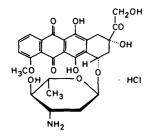


580.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	476 nm 233 nm		484 nm 232 nm	Decom- position observed
E ^{1%} 1cm	206 659		194 614	
3	12000 38200		11200 35600	

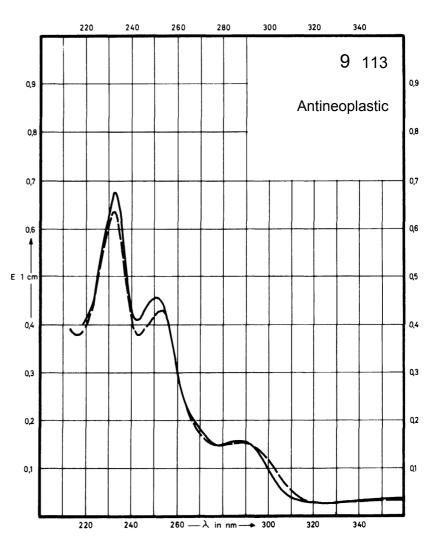

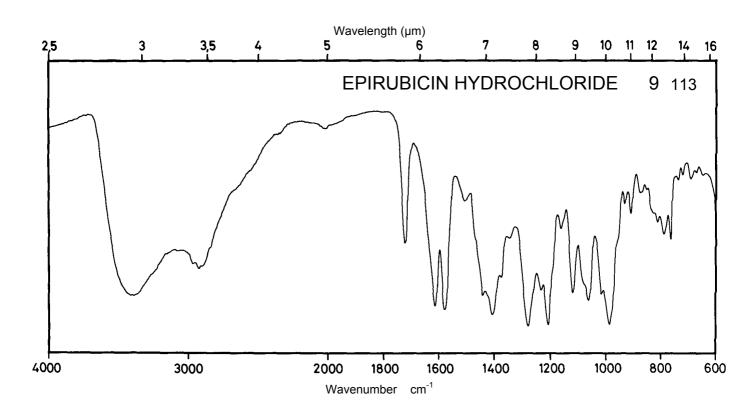



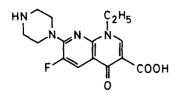



320.3

Concentration


0.84 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	344 nm 273 nm		339 nm 268 nm	346 nm 266 nm
E ^{1%} 1cm	485 1072		514 1183	544 926
ε	15500 34400		16450 37900	17400 29650

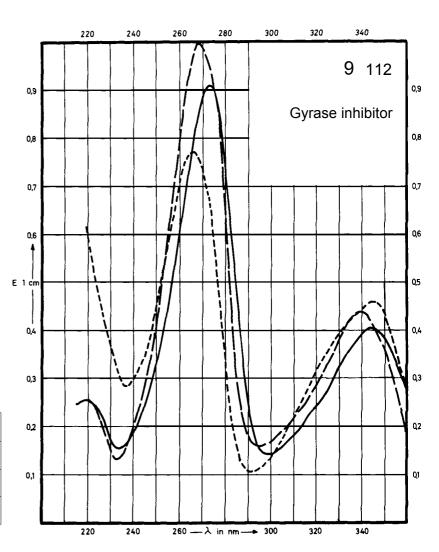

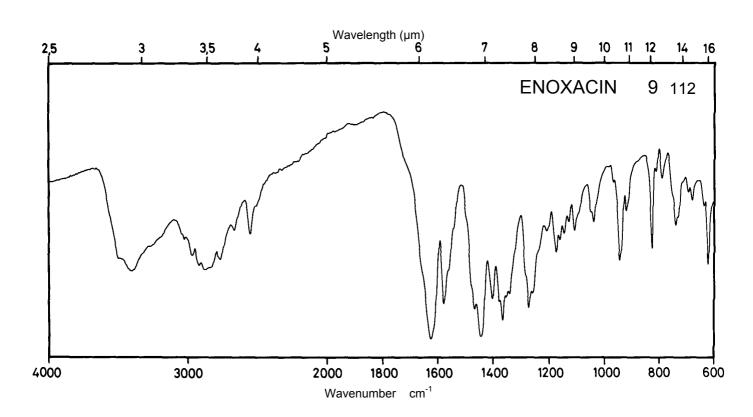




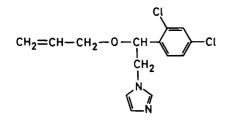


 M_{r} 297.2


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				•••••
Maximum of absorption	281 nm 273 nm			
E ^{1%} 1cm	9.5 11.4			
ε	280 340			

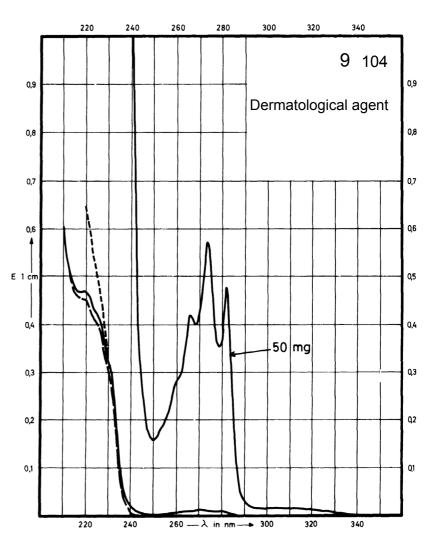

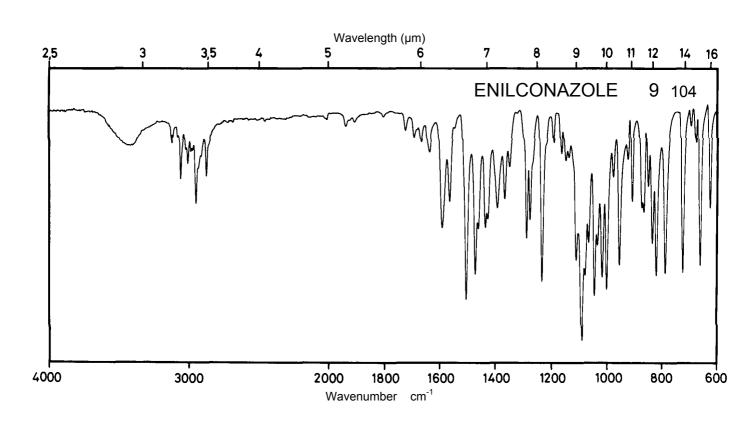




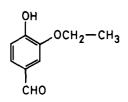


 M_{r} 166.2


Concentration 0.2 mg / 100 ml 1.0 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	308 nm 278 nm 230 nm		308 nm 279 nm 230 nm	347 nm 248 nm
E ^{1%} 1cm	580 625 945		545 645 905	1540 573
ε	9600 10400 15700		9000 10700 15100	25600 9500

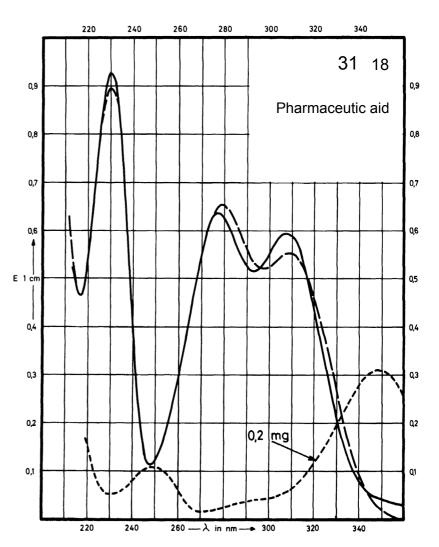

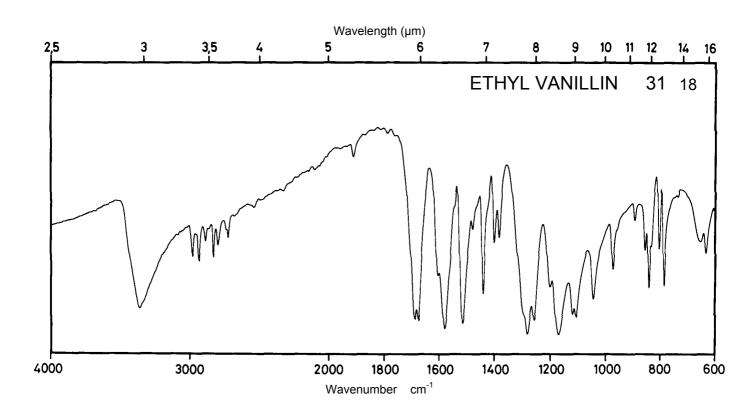


Name ETHYL BISCOUMACETATE

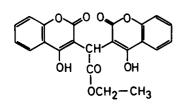


 M_{r}


408.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm 288 nm	303 nm 276 nm		311 nm
E 1% 1cm	467 424	398 449		622
ε	19090 17310	16240 18330		25420

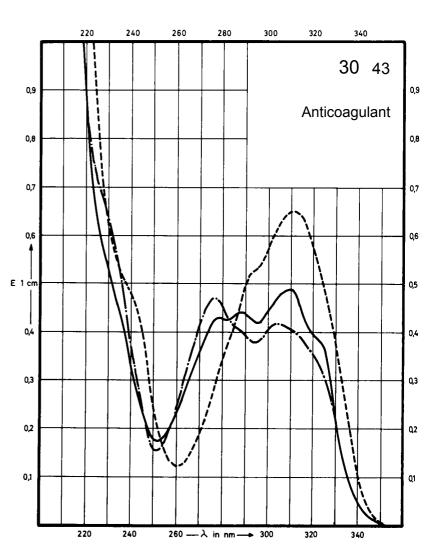

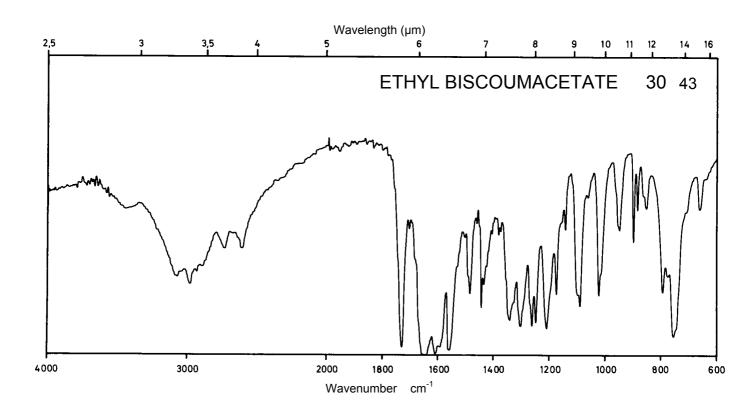





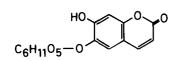


340.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	334 nm 223 nm		334 nm 223 nm	Decom- position observed
E ^{1%} 1cm	350 398		366 429	
ε	11910 13540		12450 14600	

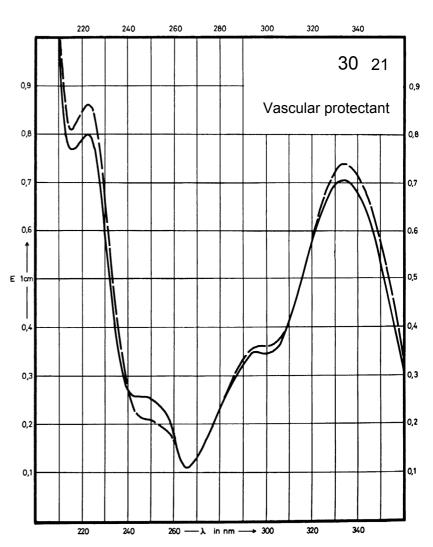

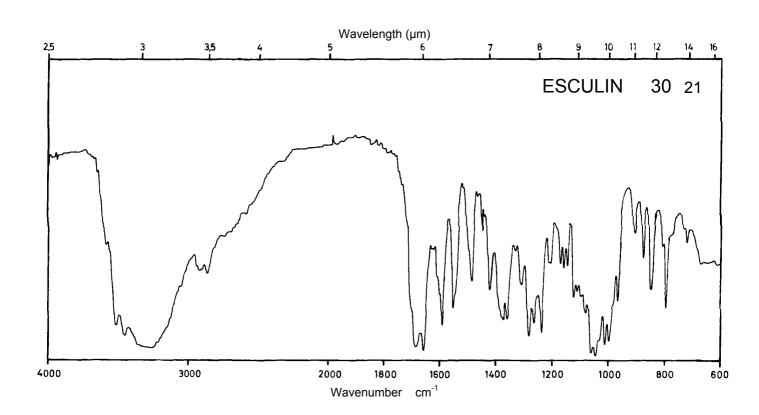

Name EPIANHYDRO-TETRACYCLINE HYDROCHLORIDE

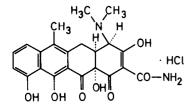


 M_{r}


462.9

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	426 nm 271 nm		426 nm 273 nm	
E ^{1%} _{1cm}	175 960		170 1015	
3	8100 44500		7900 47000	

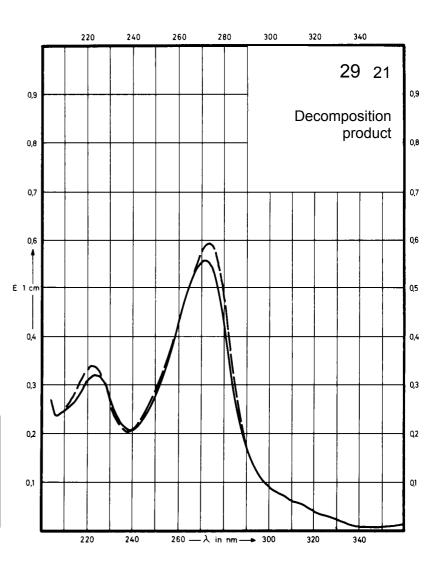

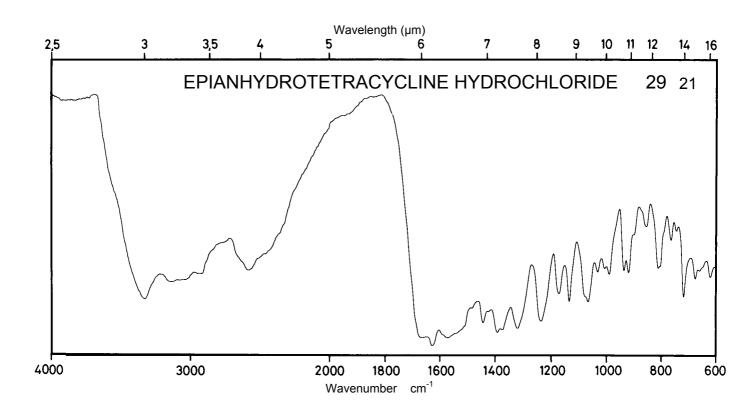


Name EPITETRACYCLINE HYDROCHLORIDE

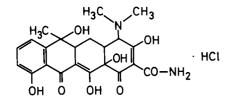


 M_r


480.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	368 nm 256 nm		359 nm 256 nm	381 nm 267 nm
E 1% 1cm	318 329		299 336	379 321
ε	15300 15800		14400 16200	18200 15400

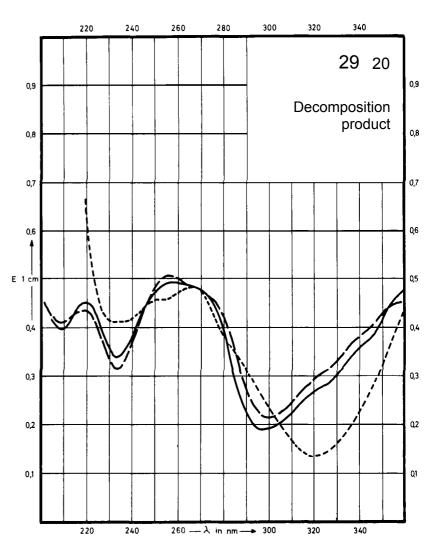

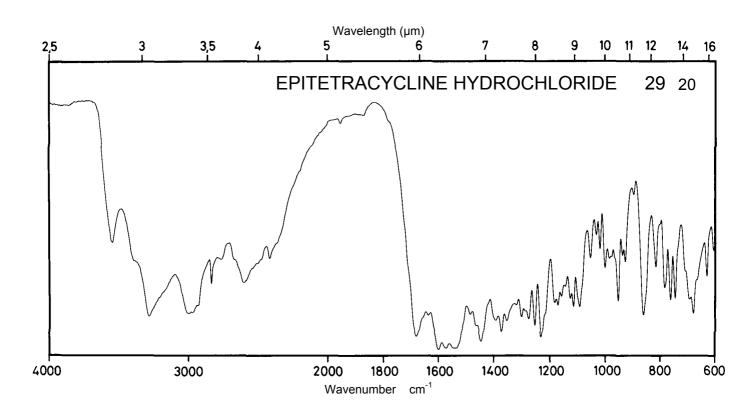

Name ETILEFRINE HYDROCHLORIDE

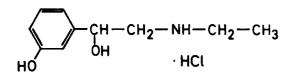



 M_{r}

217.7


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm		272 nm	290 nm 237 nm
E ^{1%} 1cm	98		84	138 409
3	2130		1830	3000 8900

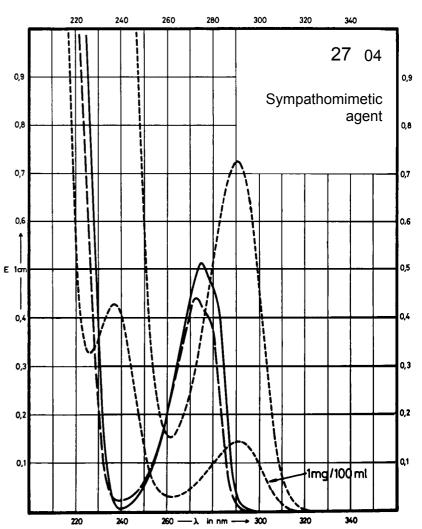

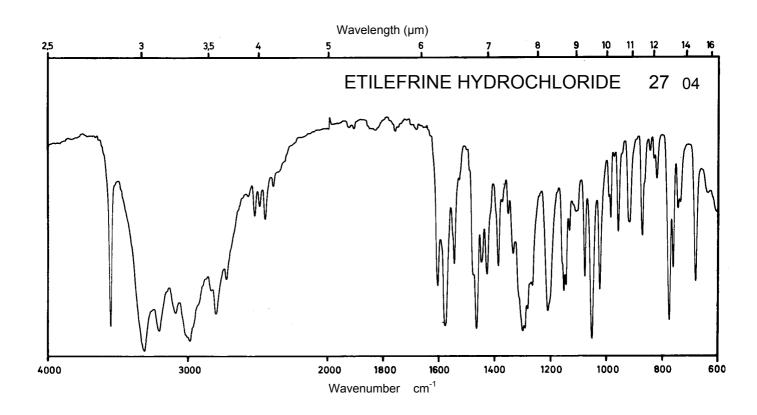





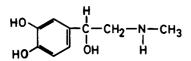


183.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm		282 nm	
E ^{1%} 1cm	162		149	
3	2970		2730	

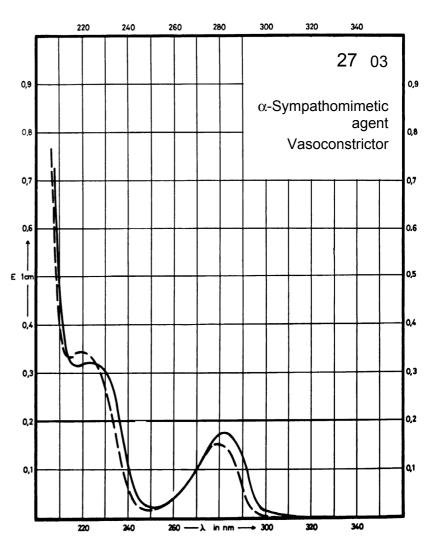

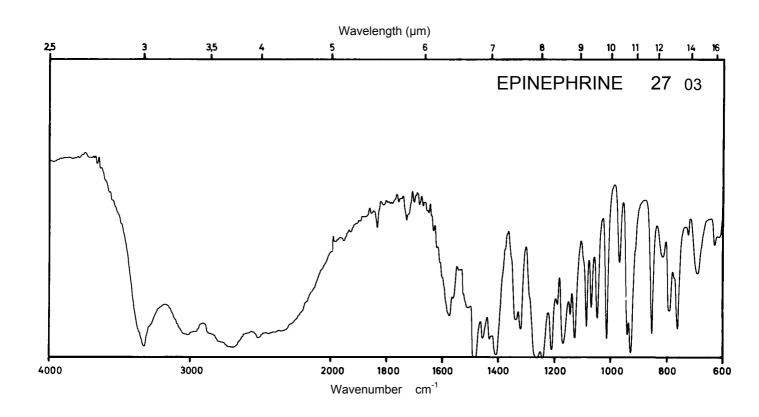




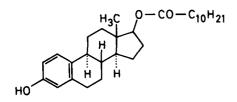


 M_{r}


440.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm			Decom- position observed
E ^{1%} 1cm	46.6			
ε	2055			

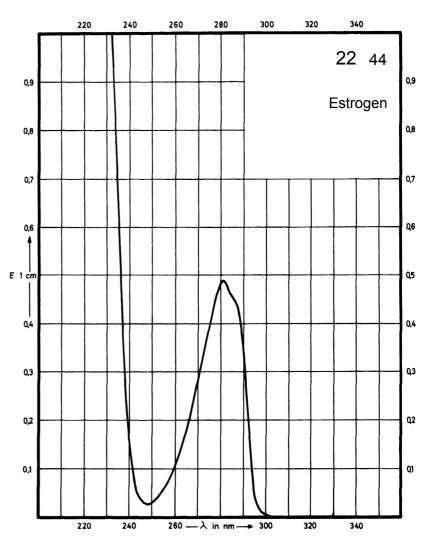

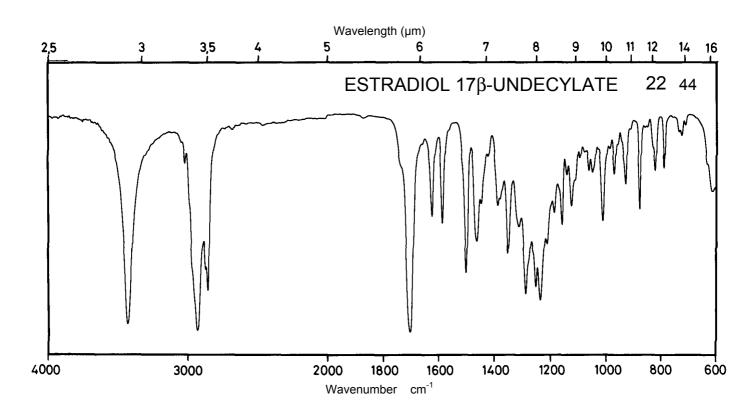




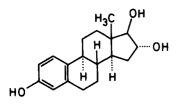


 M_r 288.4


Concentration 2 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		278 nm	296 nm 237 nm
E ^{1%} 1cm	73		66	97 319
3	2090		1920	2790 9200

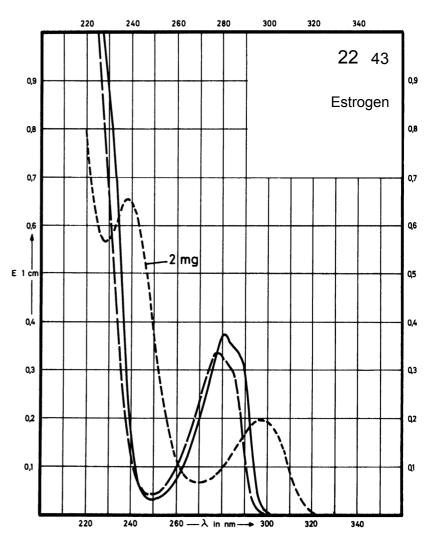

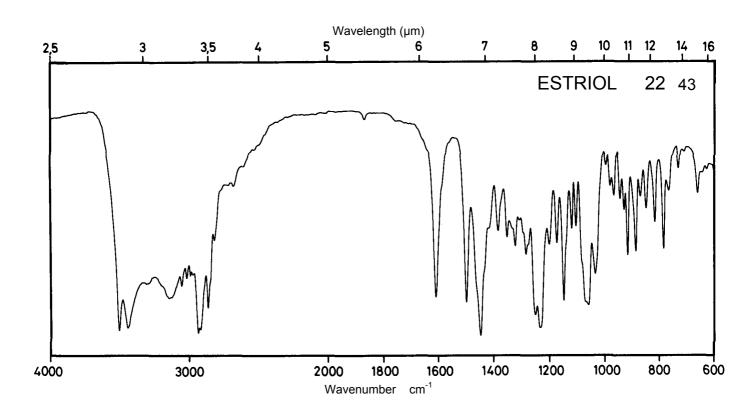





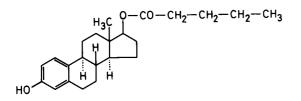


356.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm			
E ^{1%} 1cm	58			
ε	2100			

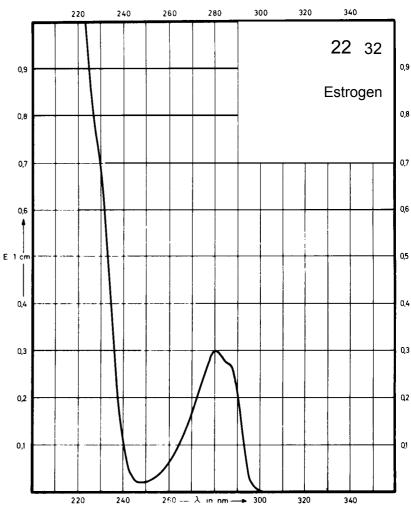

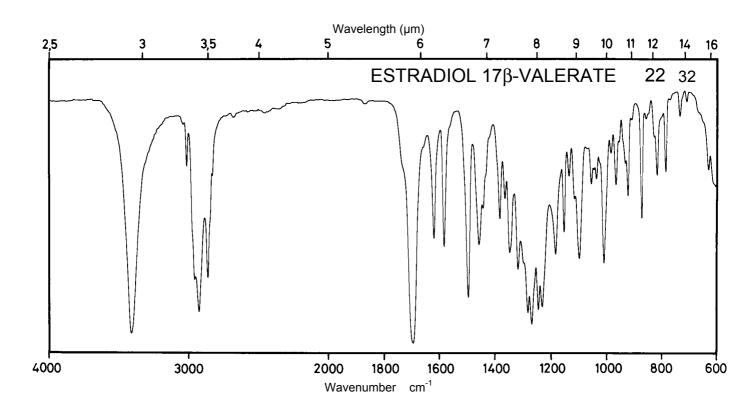





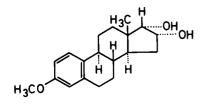


 M_r 302.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 279 nm			
E ^{1%} 1cm	62 66			
ε	1870 2000			

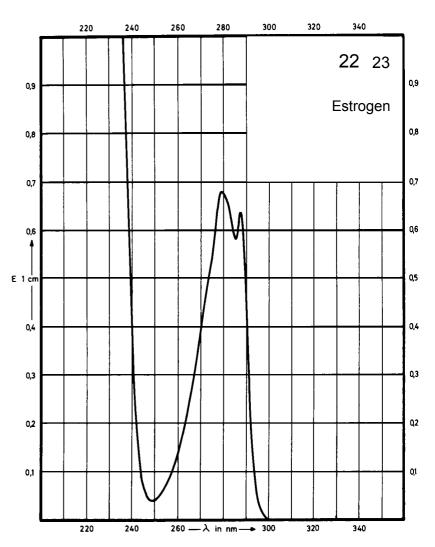

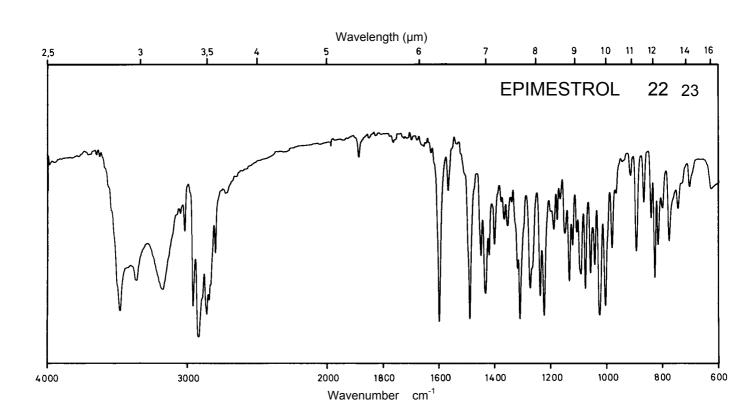





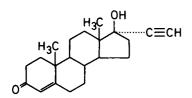


312.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm			
E ^{1%} 1cm	536			
ε	16740			

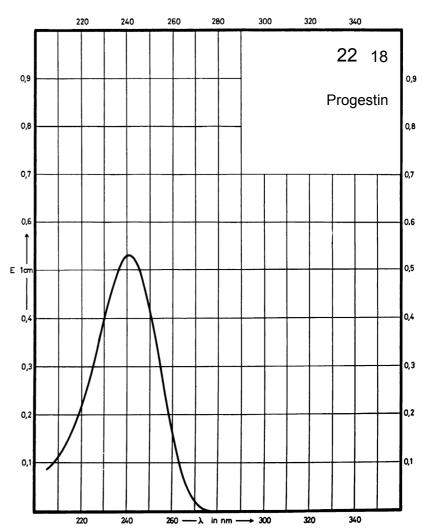

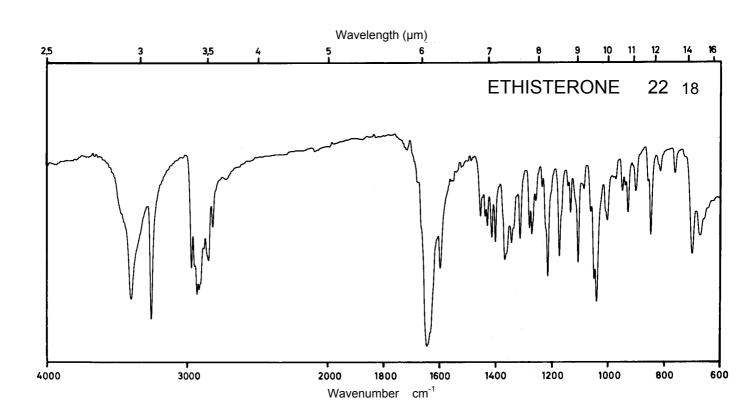





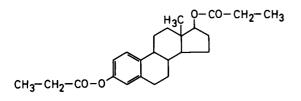


384.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 268 nm			
E ^{1%} 1cm	20 21			
ε	770 810			

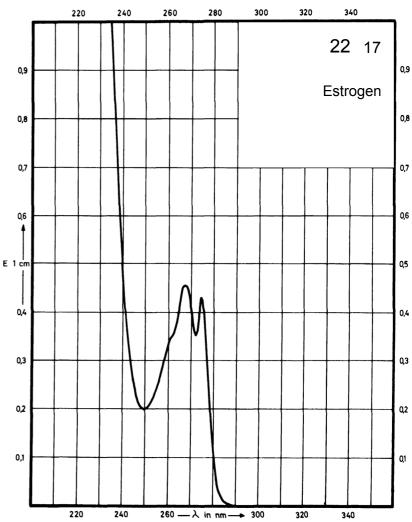

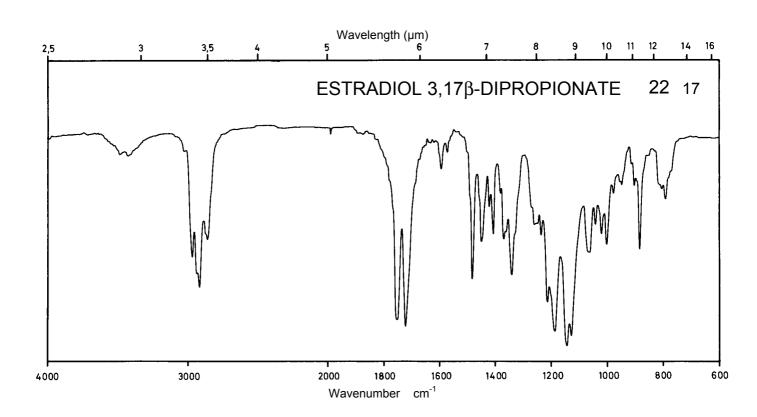




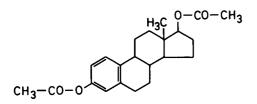


 M_{r}


356.5

F				
Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 268 nm			
E ^{1%} 1cm	20 20			
ε	710 710			

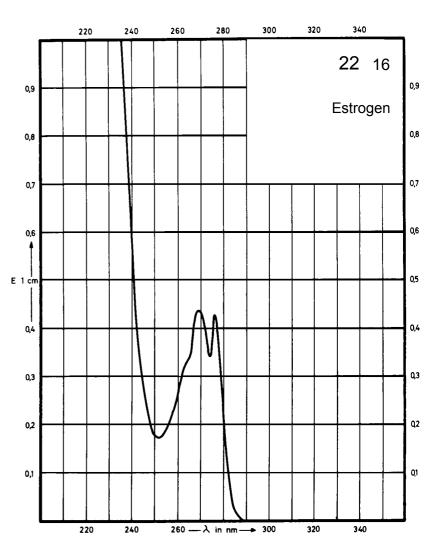

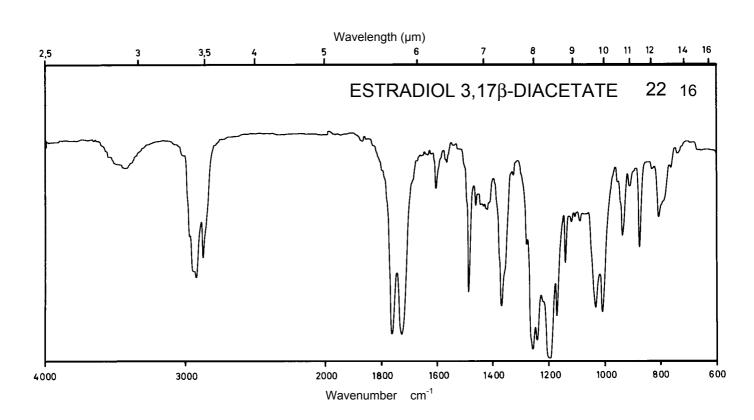

Name ESTRADIOL 3-BENZOATE

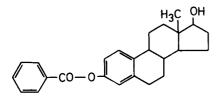



 M_{r}

376.5


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	230 nm			
E ^{1%} 1cm	510			
ε	19200			

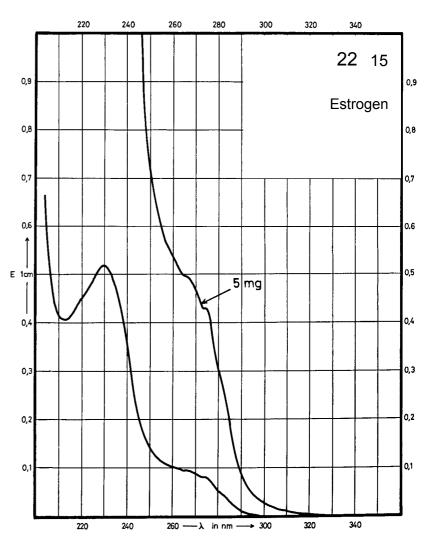

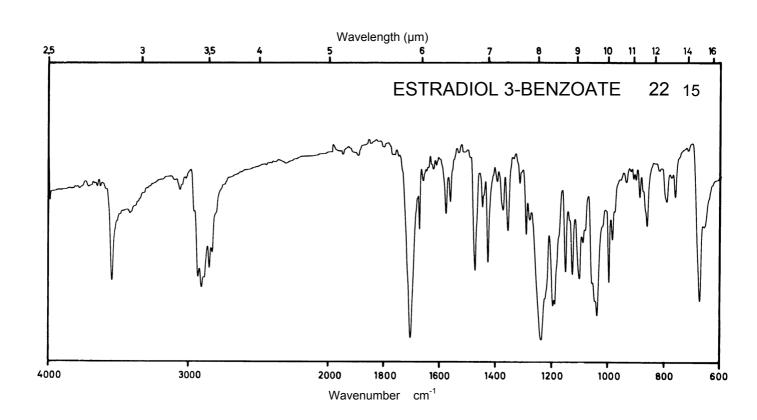


Name ETHINYLESTRADIOL

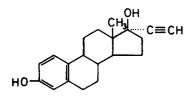


 M_{r}


296.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		278 nm	297 nm 238 nm
E ^{1%} 1cm	69		65	92 314
ε	2050		1930	2730 9310

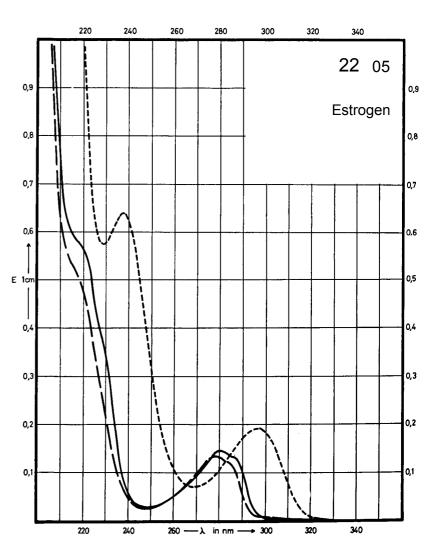

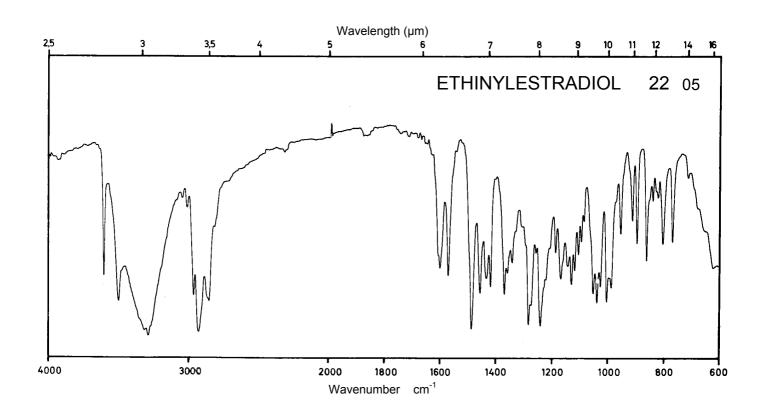




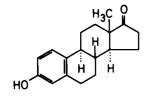


 M_r 270.4


Concentration 8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	279 nm			
E ^{1%} 1cm	78			
ε	2110			

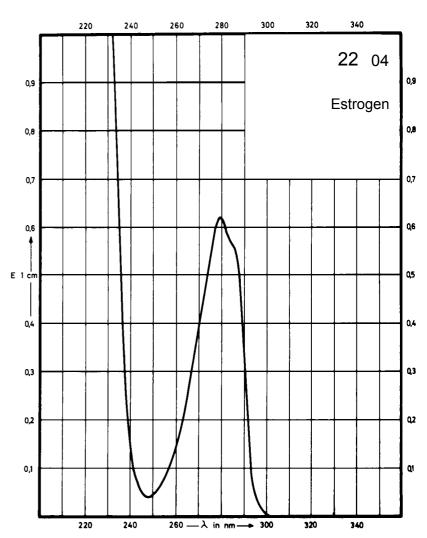

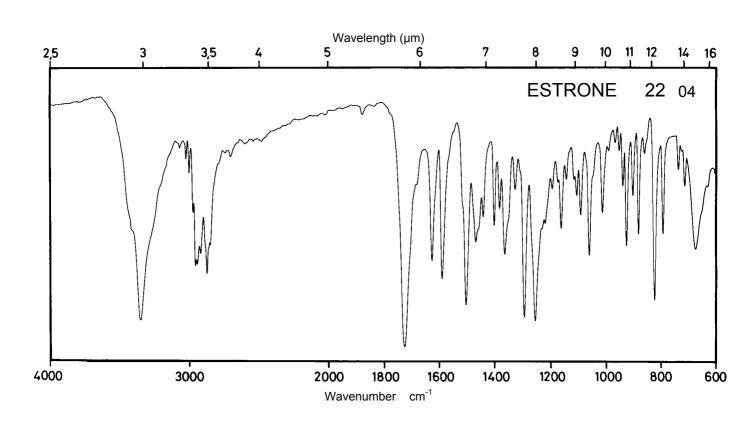

Name ESTRADIOL 17β-CIPIONATE

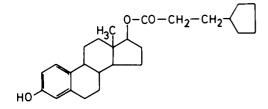



 M_{r}

396.6


Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm			
E ^{1%} 1cm	53			
ε	2100			

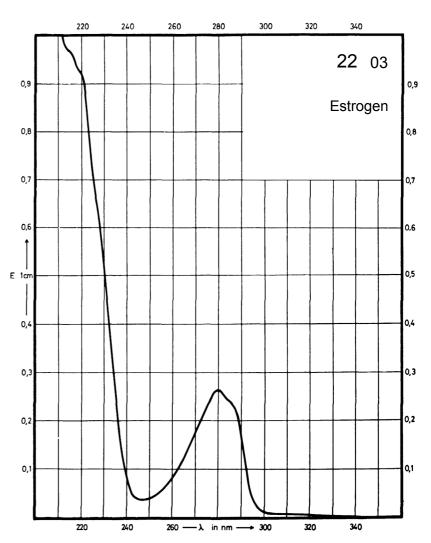

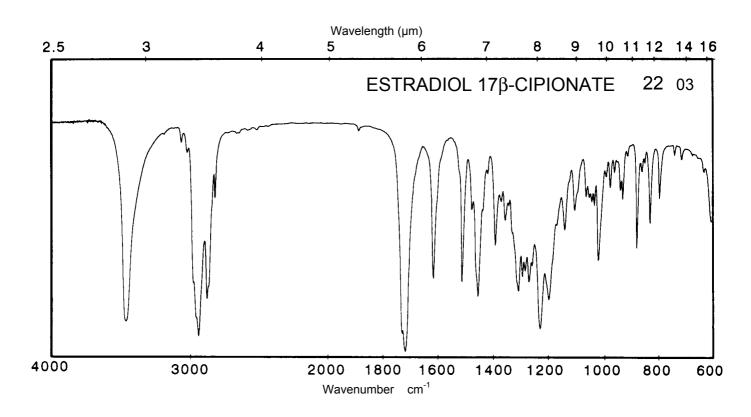



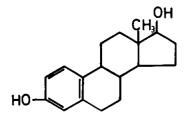



272.4

Concentration 5


5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm			
E ^{1%} 1cm	75			
ε	2040			

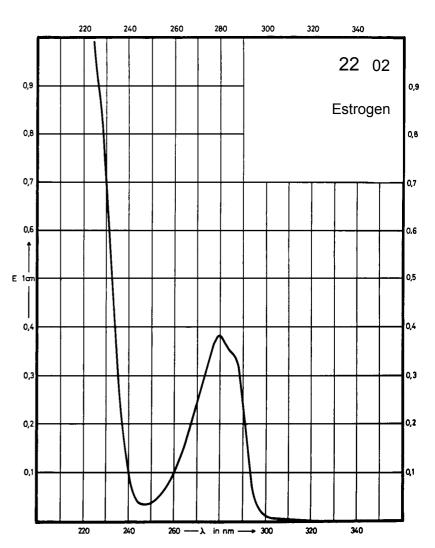

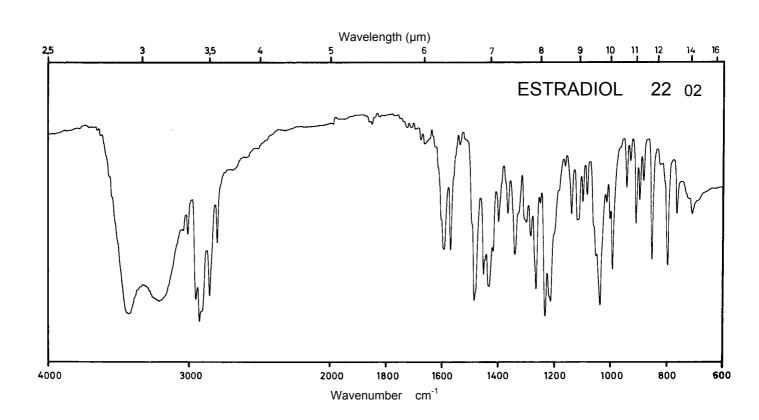




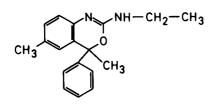


 M_r 300.8


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm	255 nm	255 nm	273 nm
E ^{1%} 1cm	449	496	538	575
ε	13500	14920	16170	17300

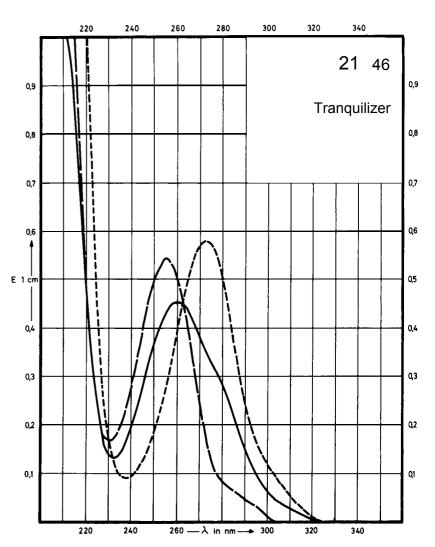

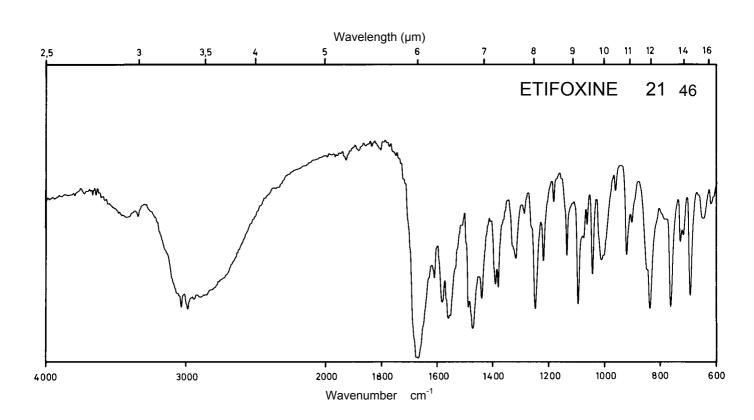



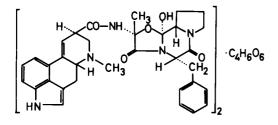



 $\,M_r\,$

1313.5


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	313 nm	318 nm	318 nm	310 nm
E ^{1%} 1cm	131	116	118	132
ε	17200	15300	15500	17400

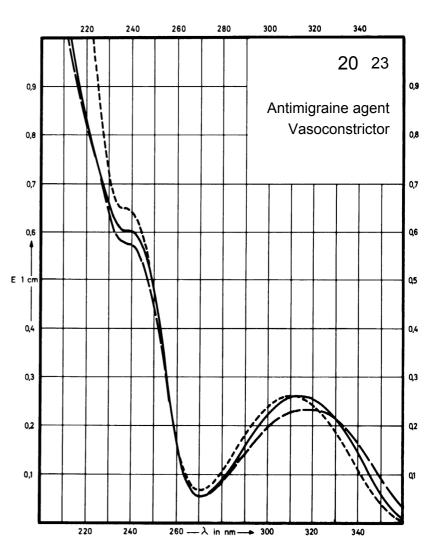

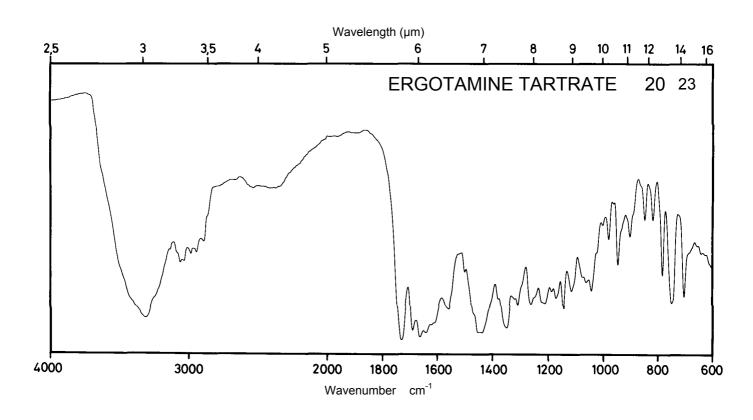



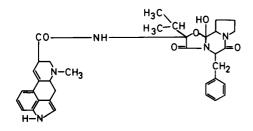



 M_r

609.7


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	312 nm		320 nm	309 nm
E ^{1%} 1cm	140		136	140
ε	8540		8290	8540

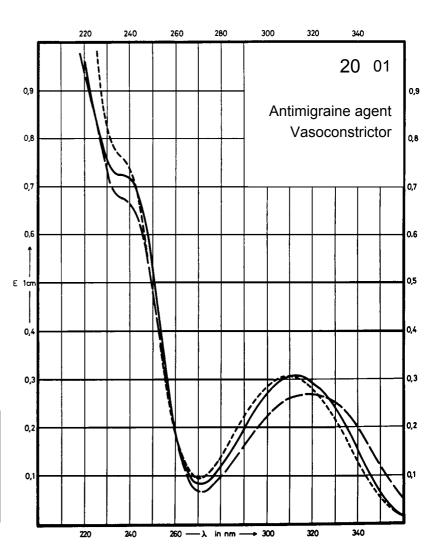

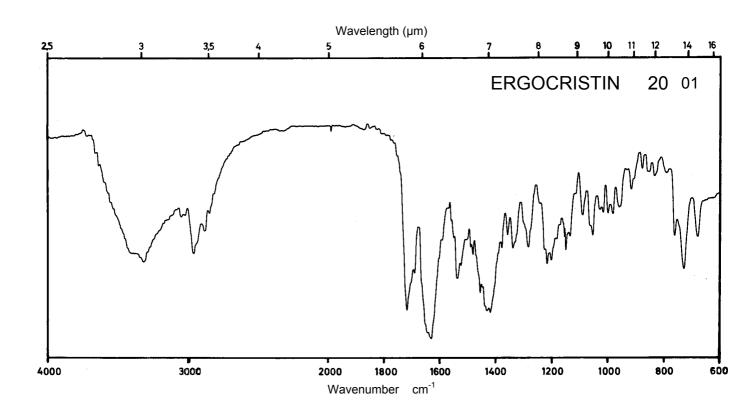




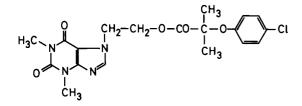


 M_r 420.9


Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		274 nm	Decom- position observed
E ^{1%} 1cm	217		212	
3	9100		8900	

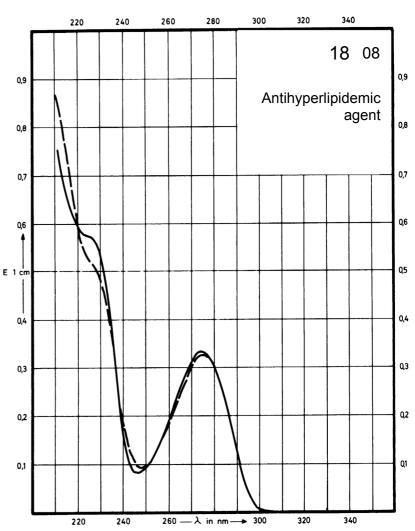

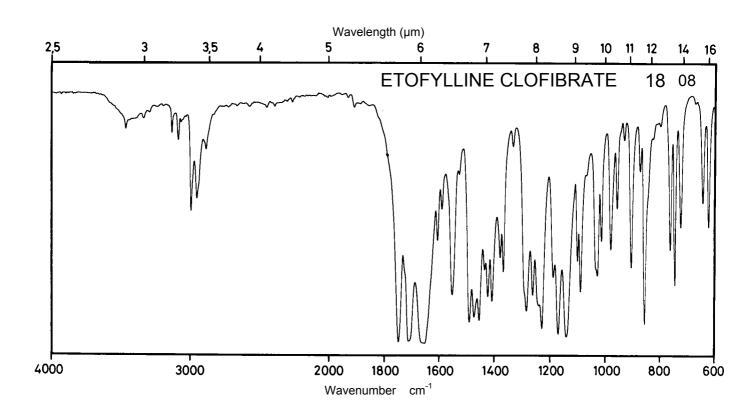



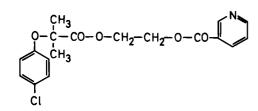

 M_r 363.8


3300 18200


ε

Concentration 1 mg / 100 ml 2 mg / 100 ml 4 mg / 100 ml


> 3300 17100

5050 13600 3150 3660

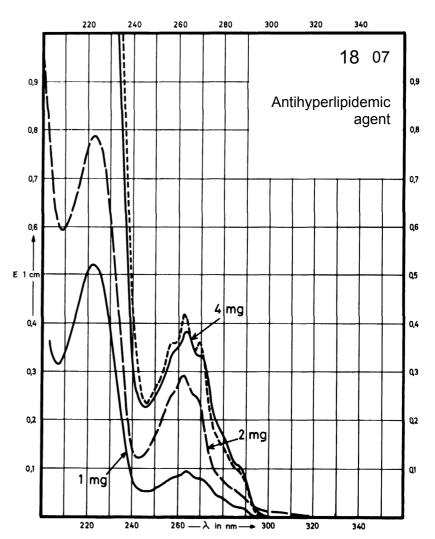

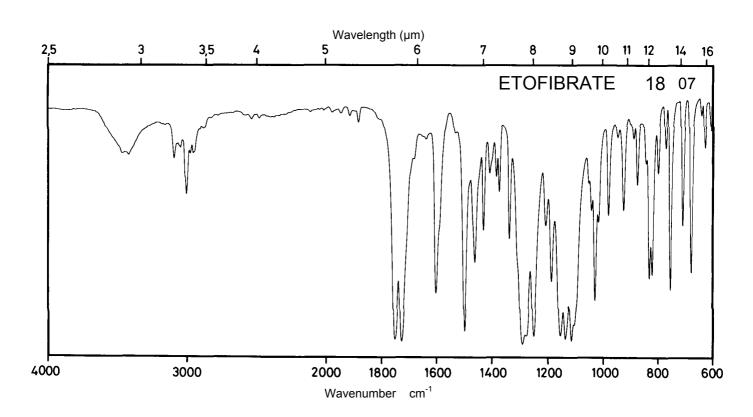




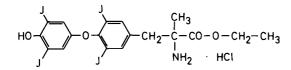


 M_r 855.4

Concentration 1 mg / 100 ml 5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	303 nm 226 nm		299 nm	327 nm
E ^{1%} 1cm	50 564		44	70
3	4260 48250		3780	5970

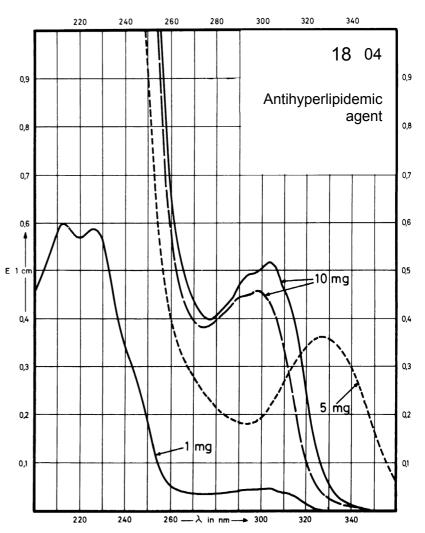

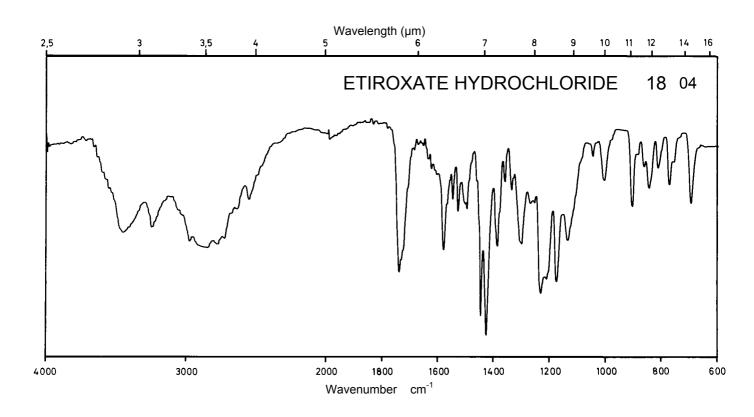


Name ETIDOCAINE HYDROCHLORIDE


 M_{r}

312.9

Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 263 nm	272 nm 263 nm	272 nm 263 nm	
E ^{1%} 1cm	10.1 13.7	12.6 15.8	12.6 15.8	
ε	315 430	396 494	396 494	

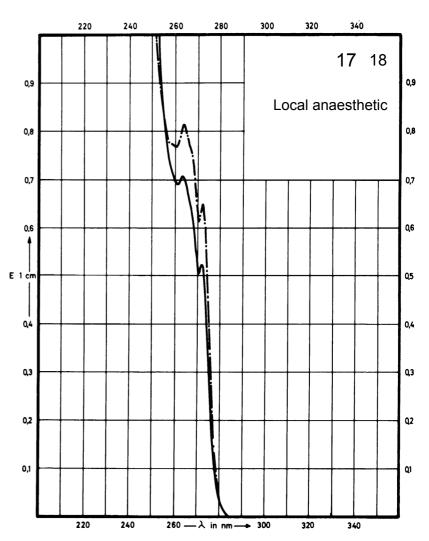

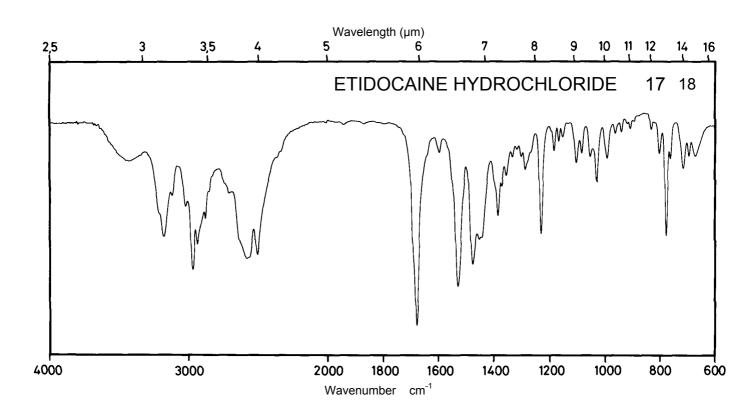





581.7


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	313 nm		318 nm	310 nm
E ^{1%} 1cm	156		135	147
3	9070		7850	8550

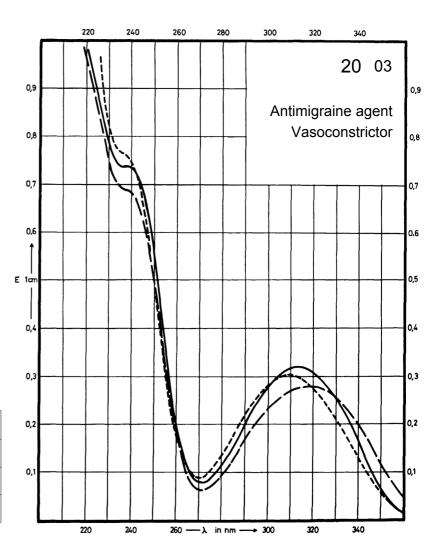

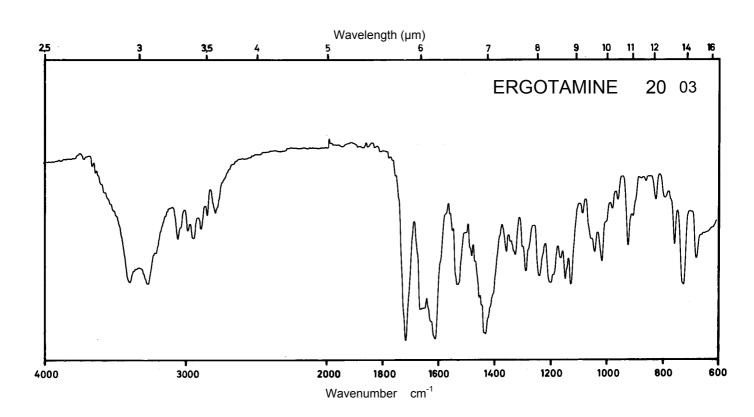




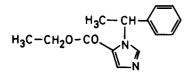


244.3

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm			
E ^{1%} 1cm	490			
ε	12000			

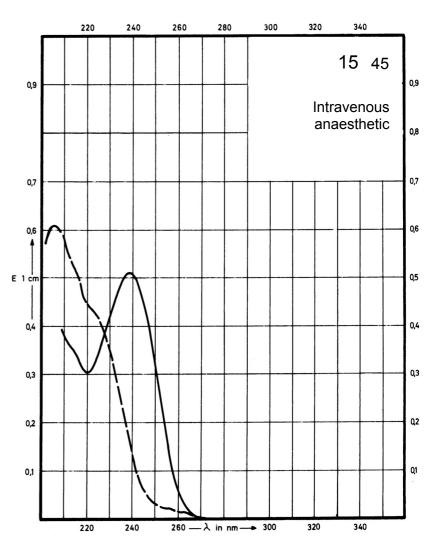

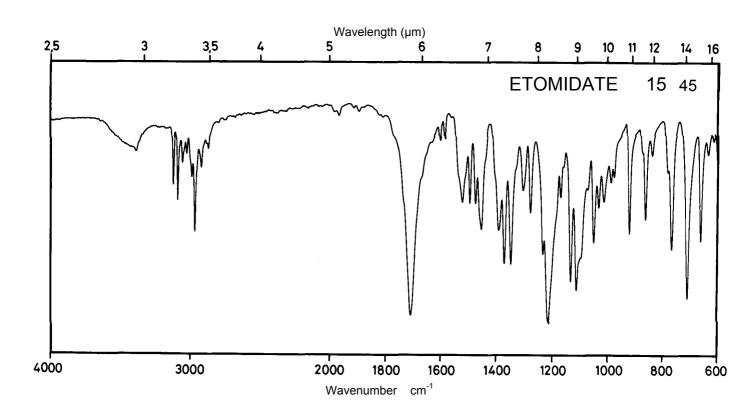


Name ETHOSUXIMIDE


 $\begin{array}{c} \text{O} \quad \text{N} \quad \text{O} \\ \text{CH}_2 - \text{CH}_3 \end{array}$


 M_r 141.2

Concentration 10 mg / 100 ml 50 mg / 100 ml


Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 249 nm 222 nm 245 nm 245 nm absorption 8.30 10.8 9.68 9.87 117 152 3 137 139

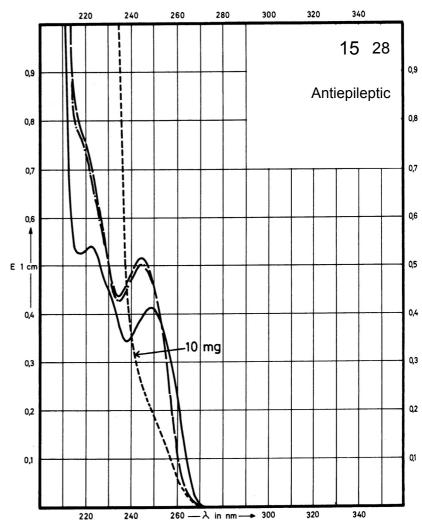


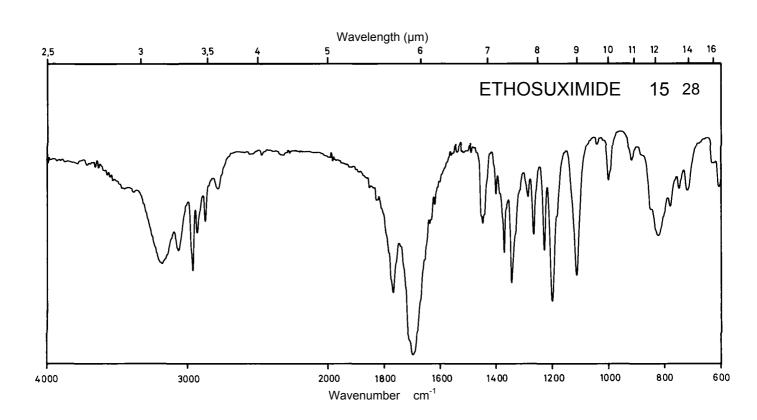

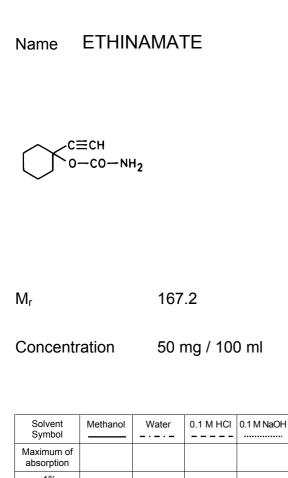


ε

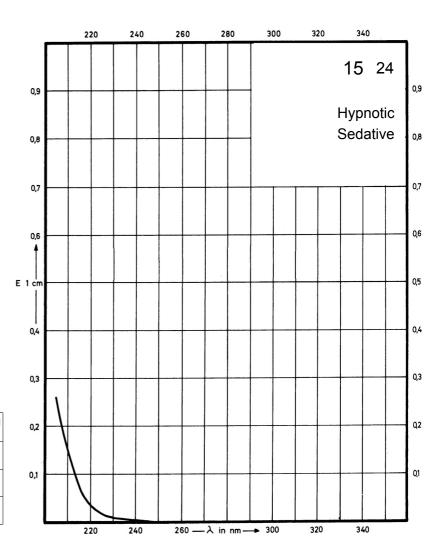

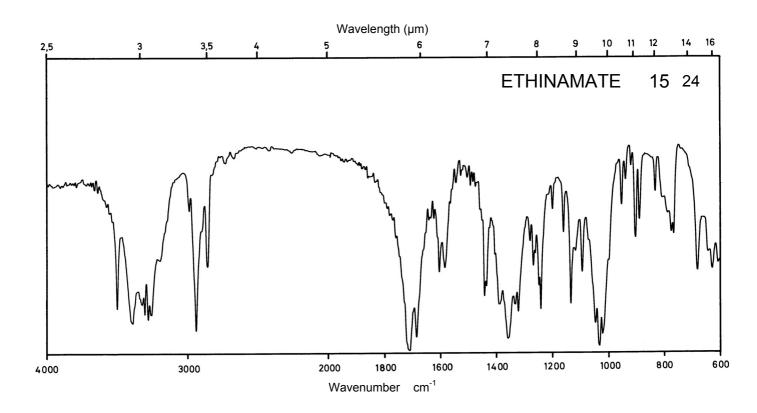




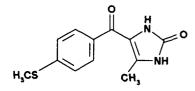


248.3


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	334 nm	333 nm	352 nm
E ^{1%} 1cm		705	719	790
ε		17500	17900	19600

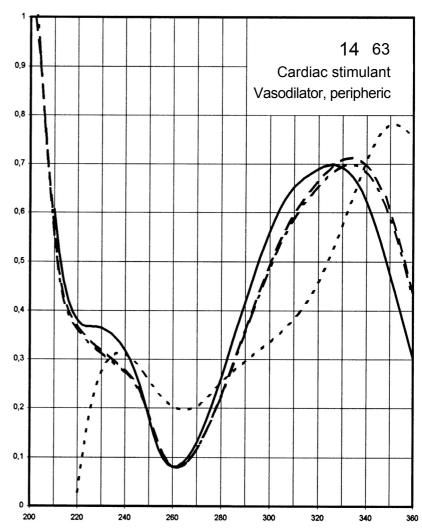

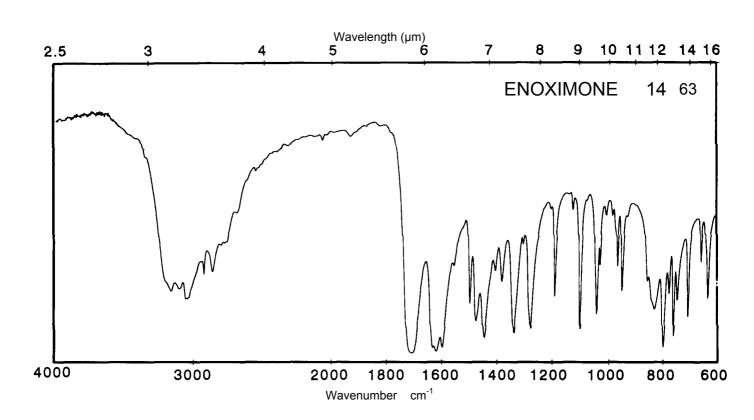



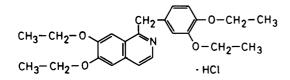

 M_r 432.0

9000 64370

ε

Concentration 0.5 mg / 100 ml 2.5 mg / 100 ml 5 mg / 100 ml


 Solvent Symbol
 Methanol
 Water
 0.1 M HCl
 0.1 M NaOH


 Maximum of absorption
 315 nm 241 nm
 313 nm 253 nm
 328 nm 281 nm 240 nm

 E 1% 1cm
 208 1490
 204 1460
 156 1450


8830 63070

6740 62640

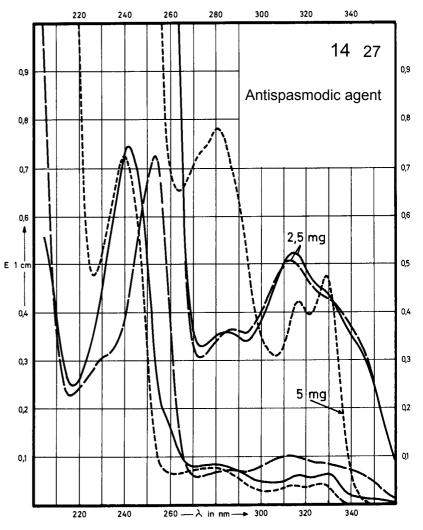

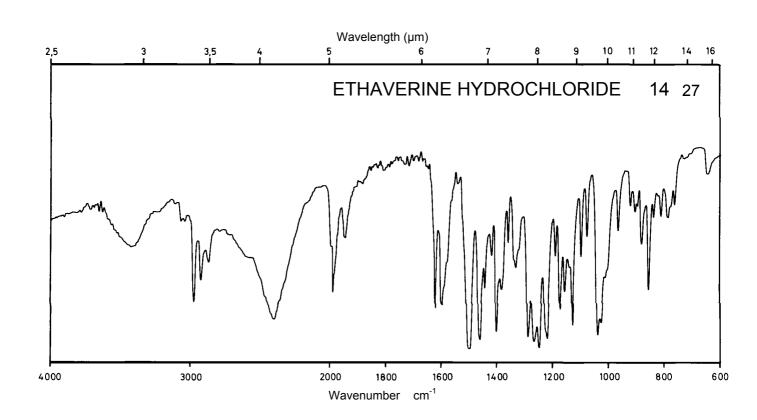




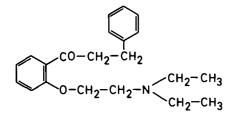


325.4

Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	302 nm 246 nm		300 nm 249 nm	307 nm 250 nm
E ^{1%} 1cm	88 213		75 193	83 184
ε	2860 6930		2440 6280	2700 5990

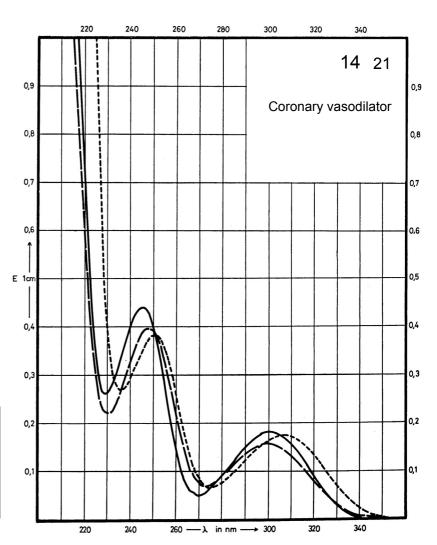

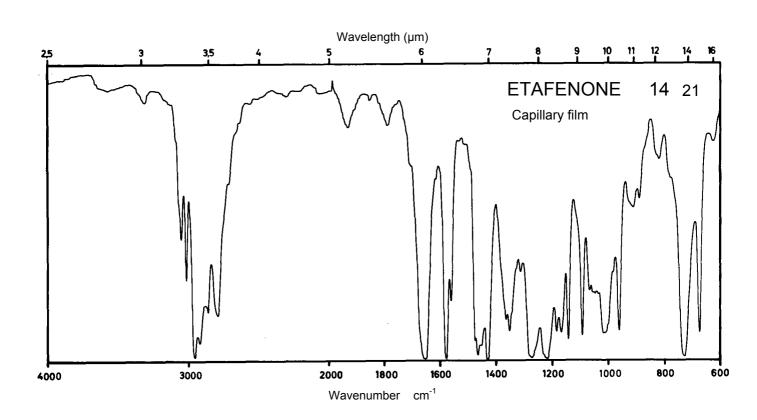

Name ETOFYLLINE


 $\begin{array}{c|c} & CH_2-CH_2OH \\ \hline \\ O = & N \\ \hline \\ CH_3 \end{array}$


 M_{r}

224.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		270 nm	272 nm
E ^{1%} 1cm	396		416	421
3	8880		9330	9440

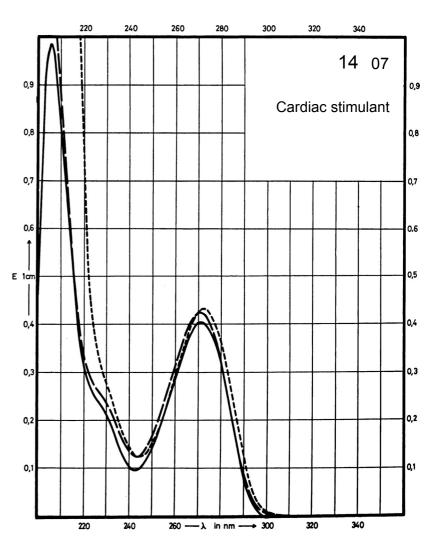

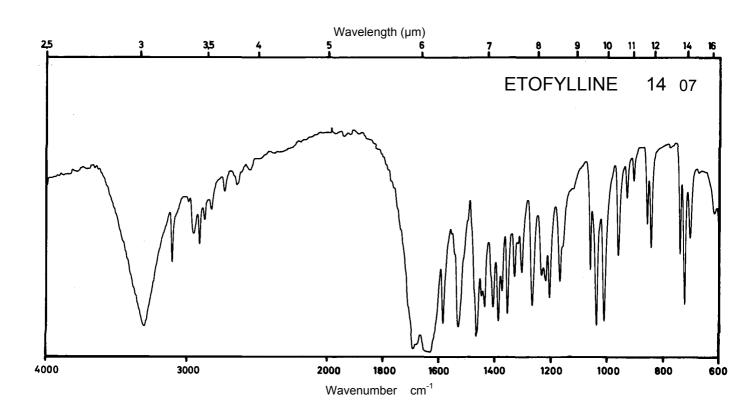




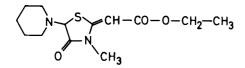


 M_r 284.4

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 245 nm	288 nm 247 nm	280 nm	Decom- position observed
E ^{1%} 1cm	741 325	752 318	740	
ε	21100 9240	21400 9040	21000	

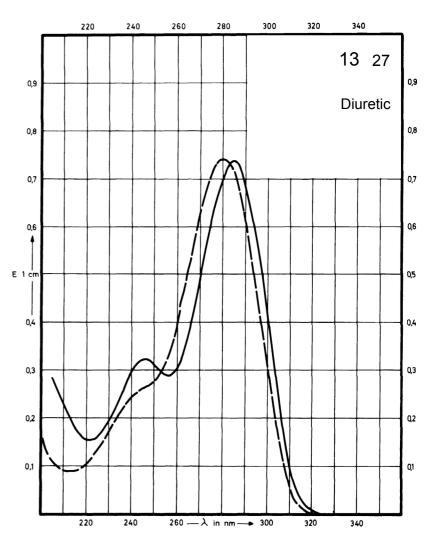

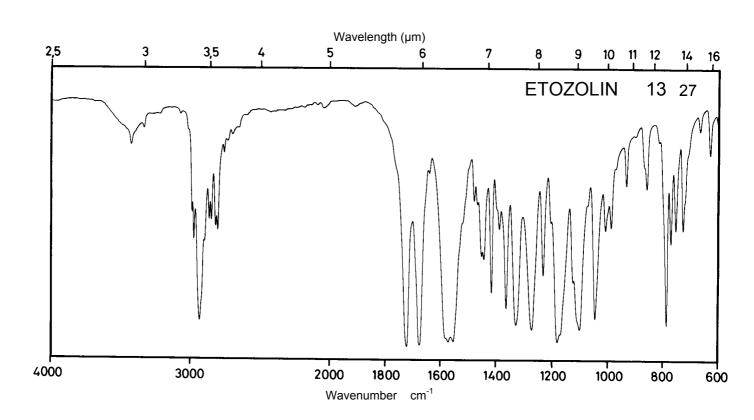


Name ETACRYNIC ACID


 M_{r}

303.1


Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm		273 nm	274 nm
E ^{1%} 1cm	103		119	163
ε	3120		3610	4940

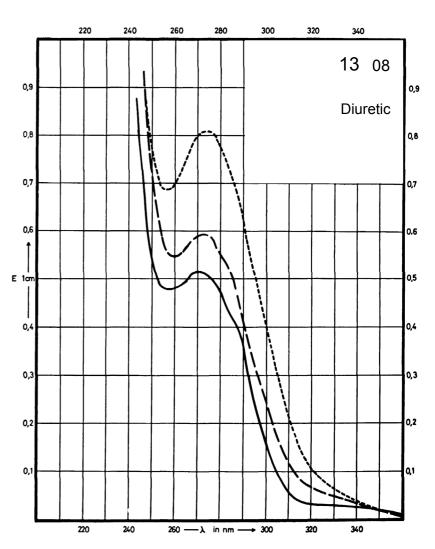

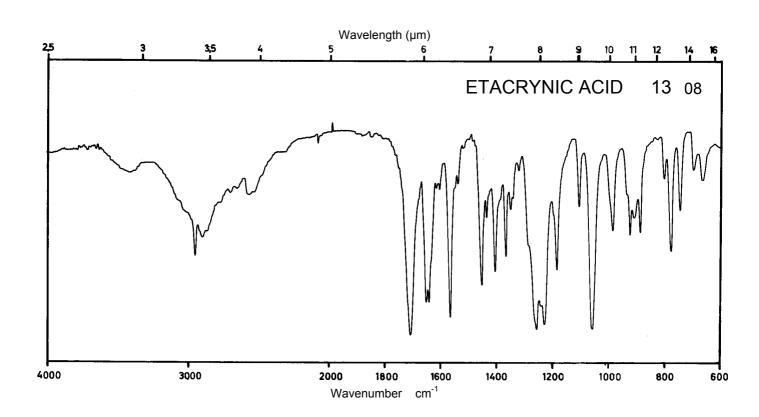



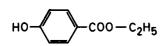



166.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	256 nm		254 nm	295 nm
E ^{1%} 1cm	1010		956	1390
з	16790		15890	23100

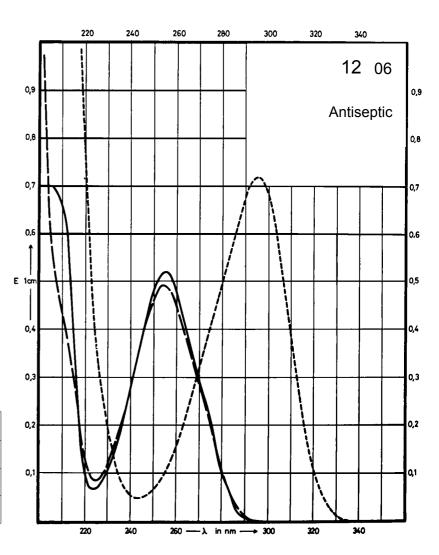

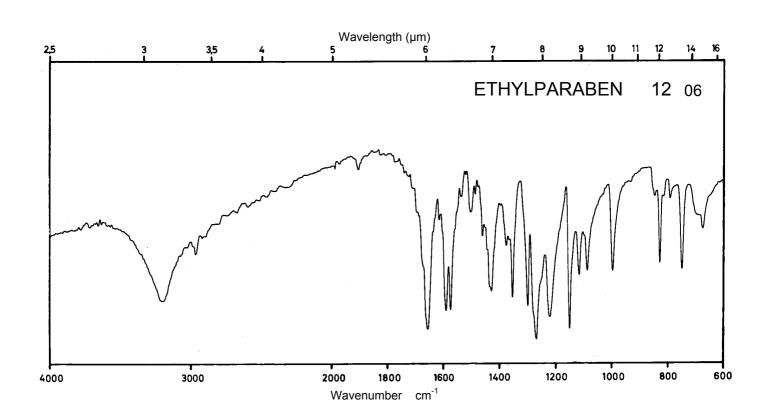



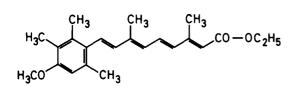



354.5

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	359 nm			
E ^{1%} 1cm	1300			
3	46100			

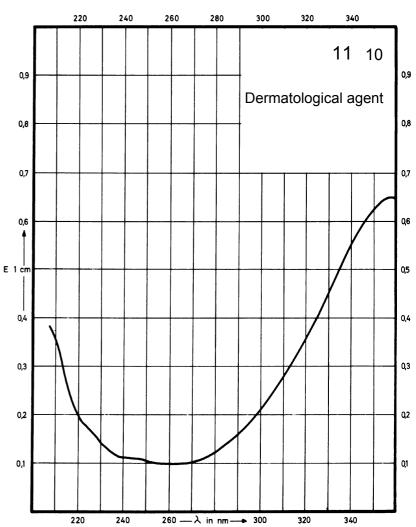

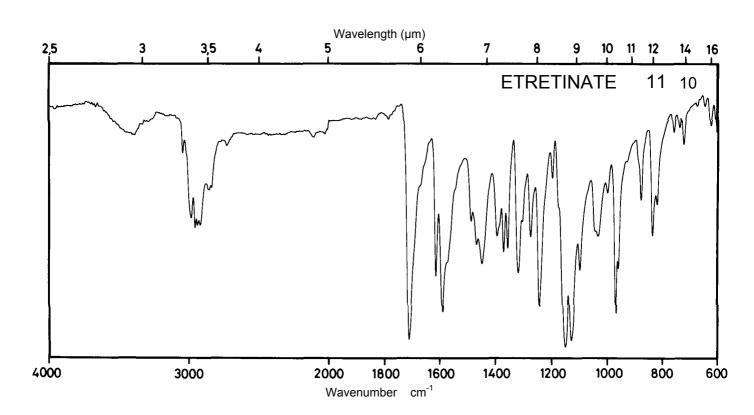




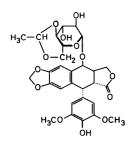



 M_r 588.6


Concentration 1 mg / 100 ml 4 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm		284 nm	287 nm 249 nm
E ^{1%} 1cm	71		67	138 212
ε	4200		3900	8100 12500

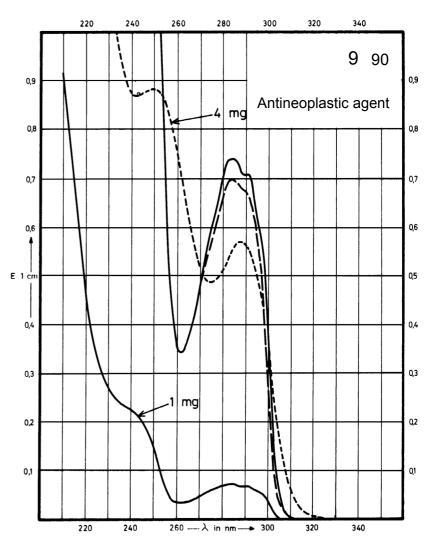


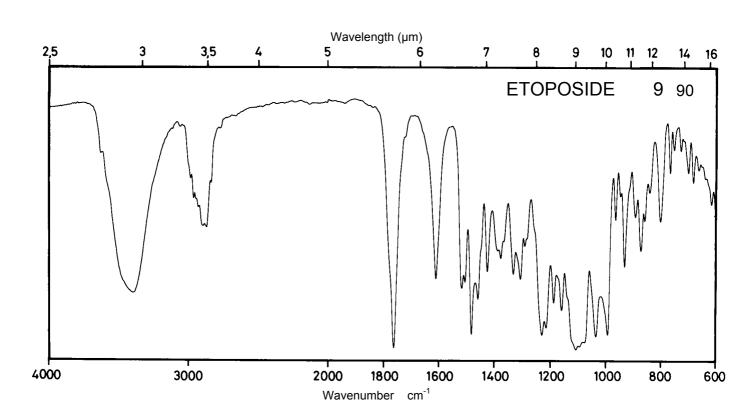

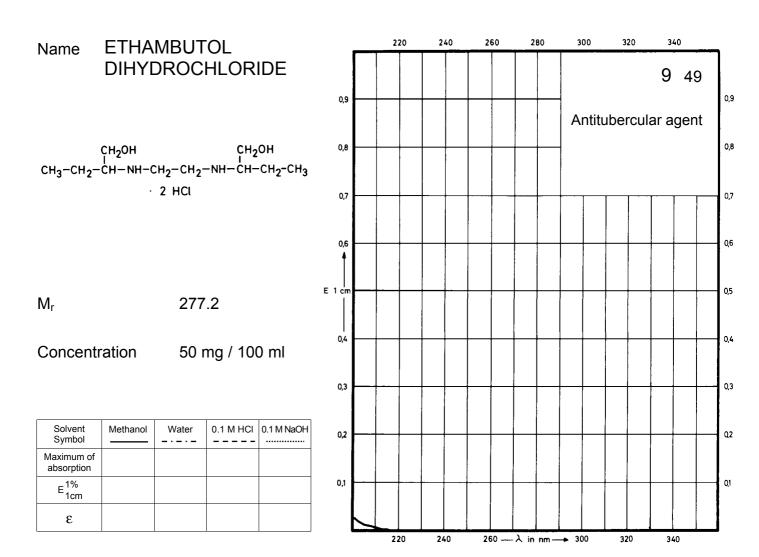

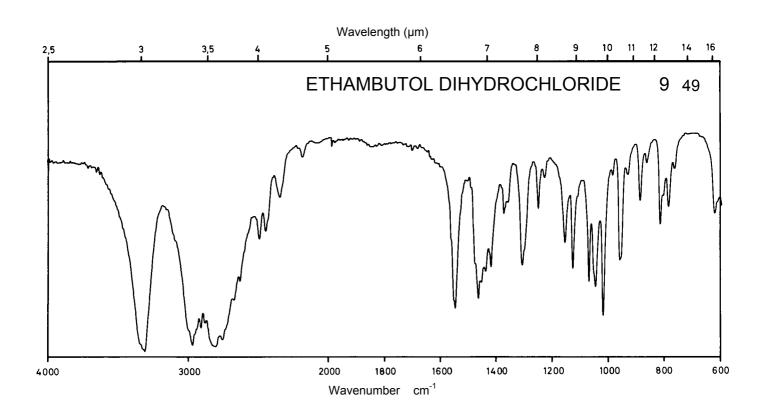




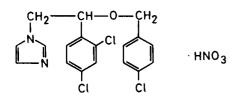


 M_r 444.7

Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm 271 nm		280 nm 271 nm	
E ^{1%} 1cm	6.5 10.8		5.9 10.4	
ε	290 480		262 462	

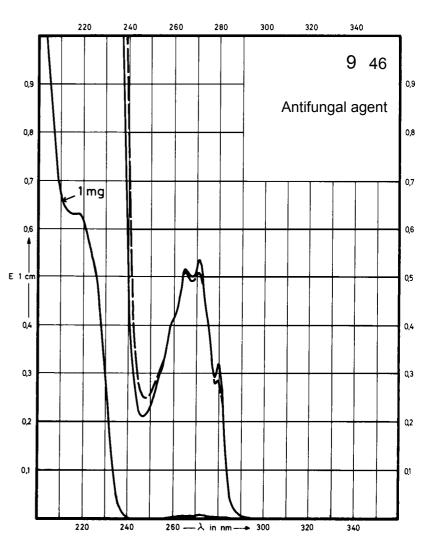

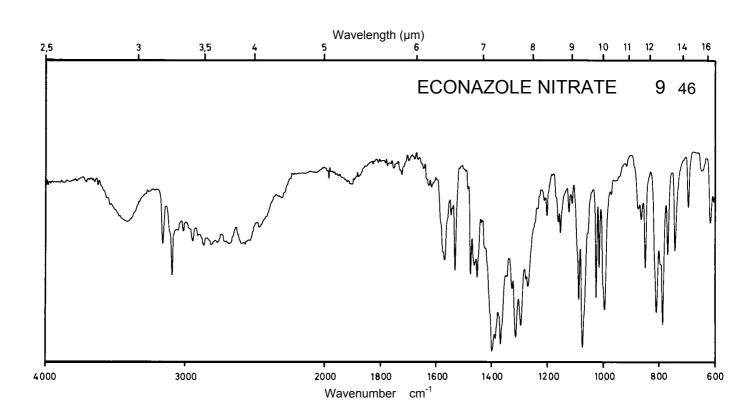







166.2

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm		277 nm 231 nm	277 nm
E ^{1%} 1cm	435		409 633	351
ε	7230		6800 10520	5840

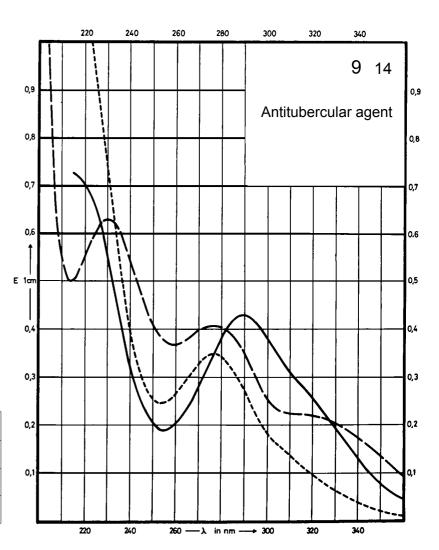

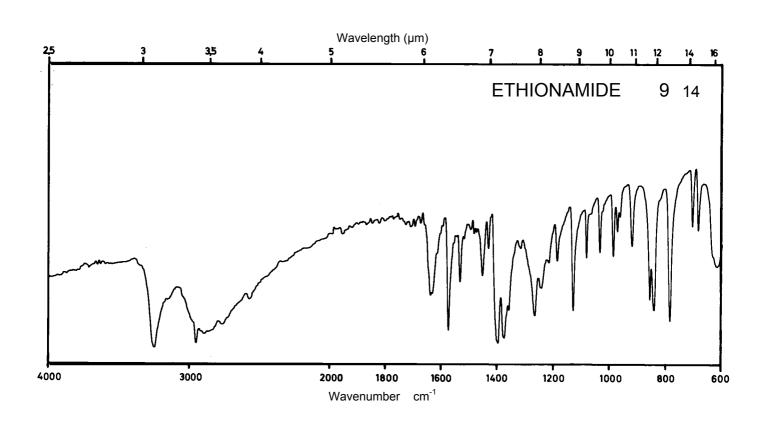



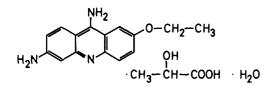



361.4

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	430 nm	422 nm	408 nm	412 nm
	373 nm	365 nm	363 nm	357 nm
	271 nm	269 nm	269 nm	274 nm
E ^{1%} 1cm	242	189	209	177
	534	412	344	239
	1760	1580	1580	1880
ε	8700	6800	7600	6400
	19300	14900	12400	8600
	63600	57100	57100	68000

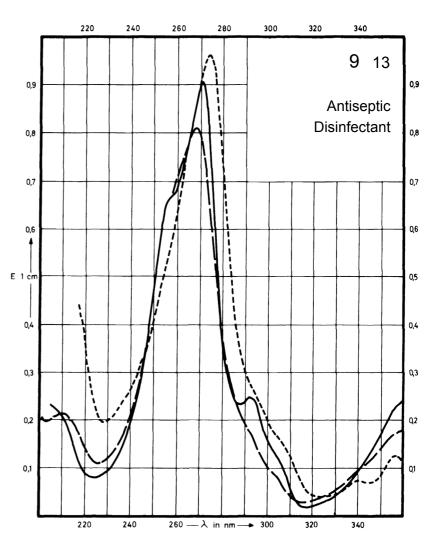

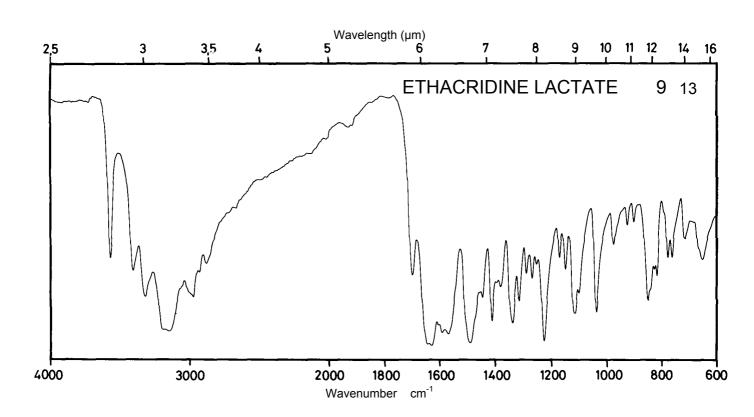

Name EMETINE DIHYDROCHLORIDE

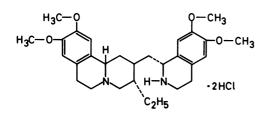


 M_{r}


553.6

Concentration


2.3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 231 nm		282 nm 230 nm	282 nm
E ^{1%} 1cm	129 289		125 286	134
з	7140 16000		6920 15830	7120

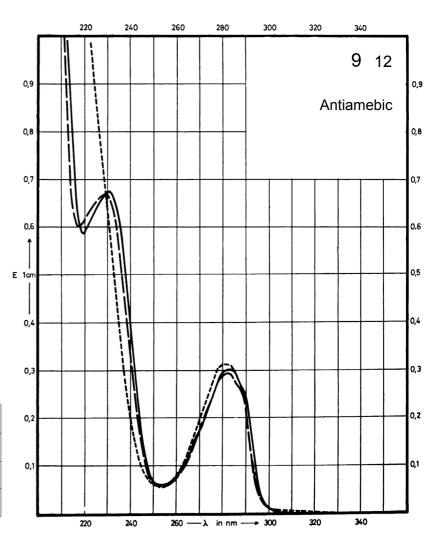

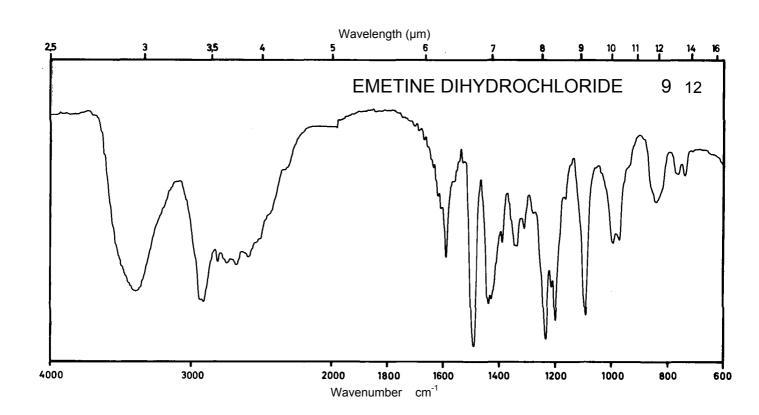




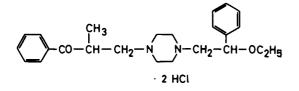


453.5


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm	249 nm	250 nm	Decom- position observed
E ^{1%} _{1cm}	267	269	282	
ε	12100	12190	12770	

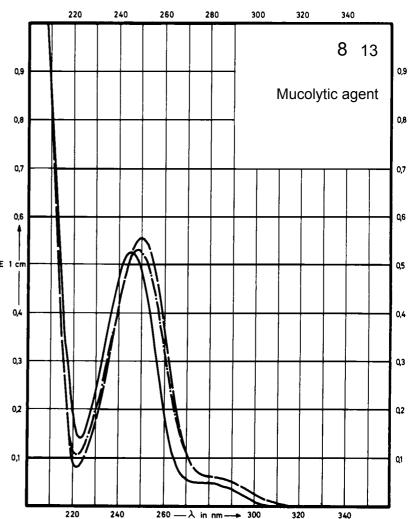

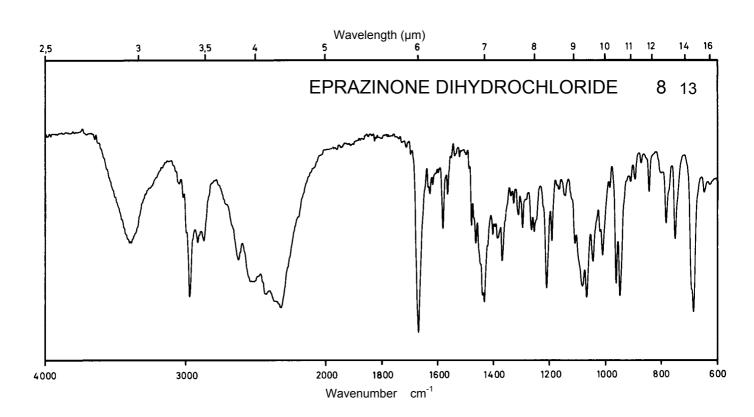

Name ETHYLMORPHINE HYDROCHLORIDE

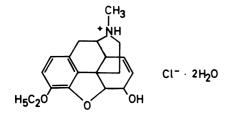



 M_{r}

385.9


Concentration 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	286 nm		284 nm	284 nm
E ^{1%} 1cm	46		44	40
ε	1780		1700	1540

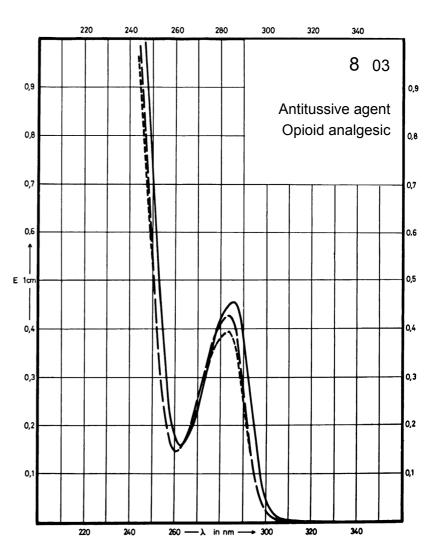

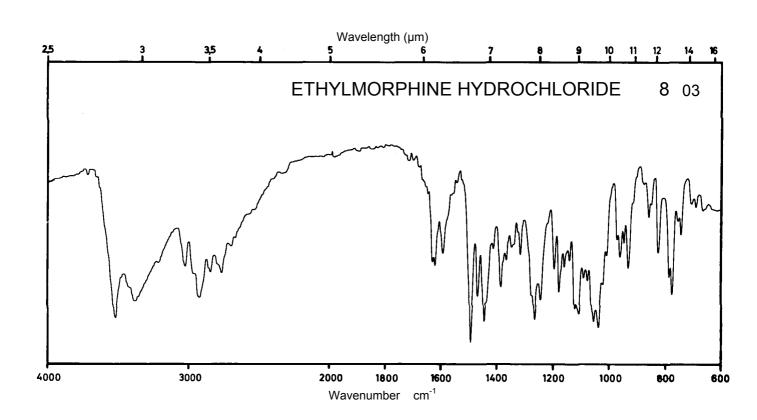



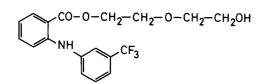



369.3

Concentration


1.1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	347 nm 286 nm	348 nm 283 nm	349 nm 284 nm	Decom- position observed
E ^{1%} 1cm	238 429	185 320	165 287	
ε	8800 15530	6820 11800	6100 10580	

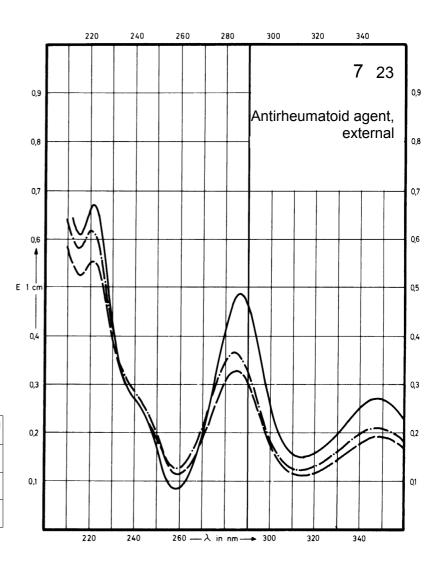

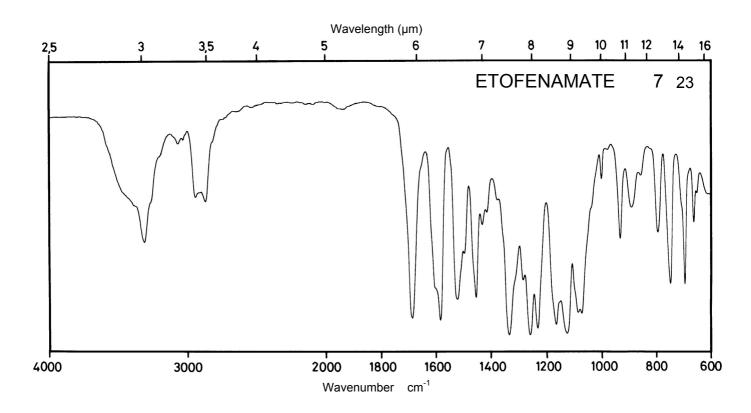




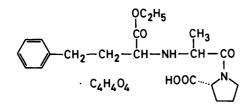


 M_r 492.5


Concentration 30 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm 264 nm 257 nm		(257 nm)	
E ^{1%} 1cm	18 22.5 30		19	
3	900 1100		930	

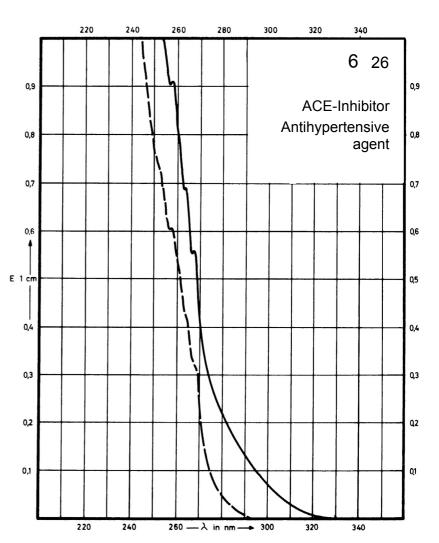

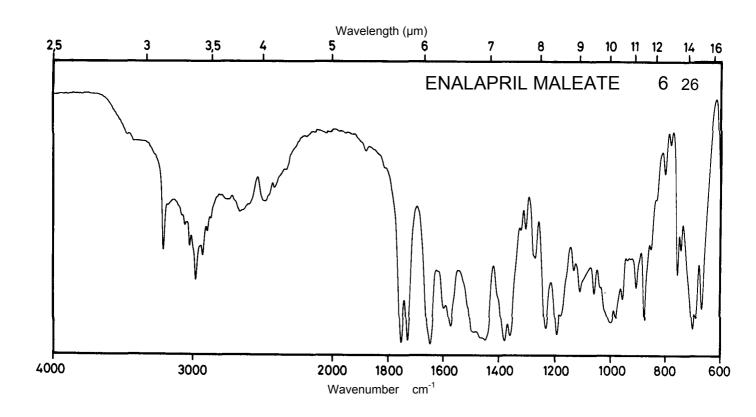




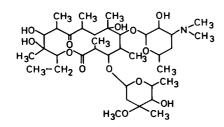


734.0


Concentration 400 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

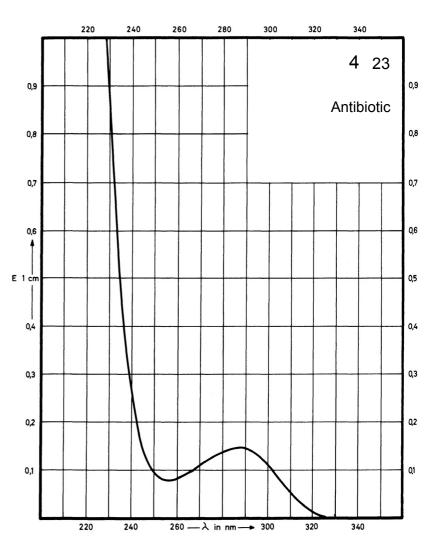

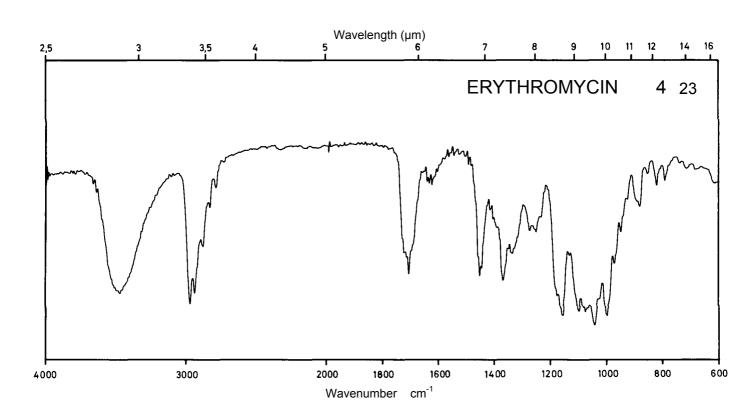



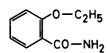



165.2

Concentration


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	290 nm 230 nm		293 nm 234 nm	291 nm
E ^{1%} 1cm	210 532		195 460	195
3	3470 8790		3220 7600	3220

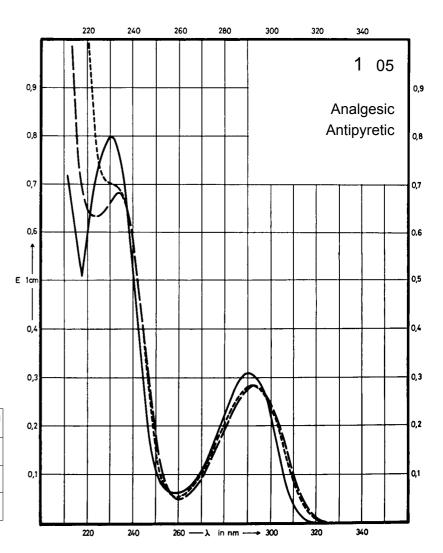

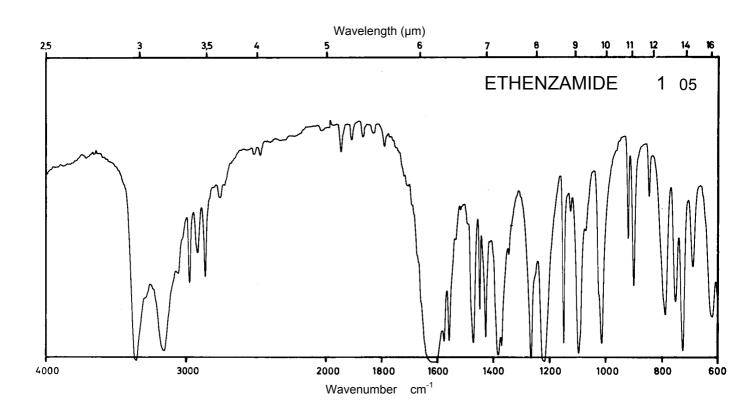



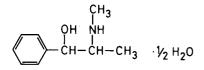



 M_r

Concentration 65 mg / 100 ml


174.2

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	264 nm		262 nm	
	258 nm		256 nm	
absorption	252 nm		250 nm	
1%	8.4		8.4	
E ^{1%} 1cm	10.9		10.9	
	8.9		8.4	
3	146		146	
	190		189	
	454	l	1.17	1

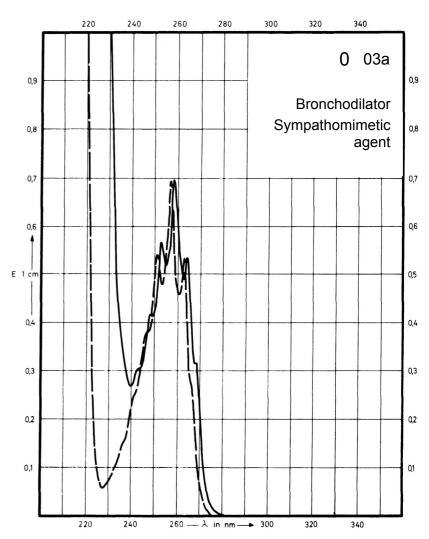

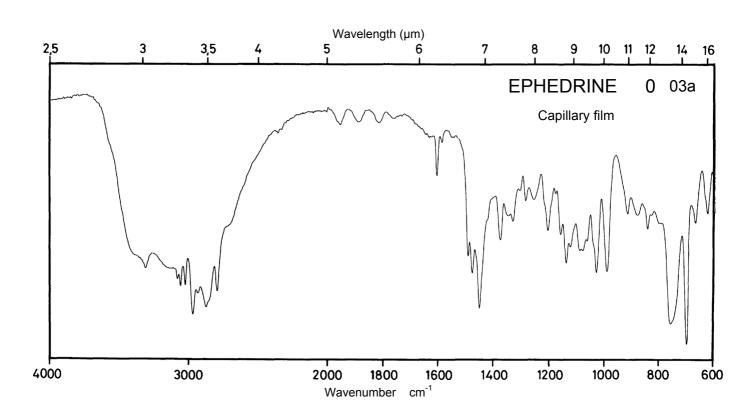




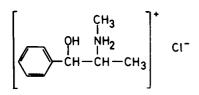


201.7


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm 256 nm 250 nm		262 nm 256 nm 250 nm	263 nm 257 nm 251 nm
E ^{1%} 1cm	8.2 10.8 8.4		7.2 9.4 7.4	7.4 9.7 7.8
ε	165 220 170		145 190 150	150 195 160

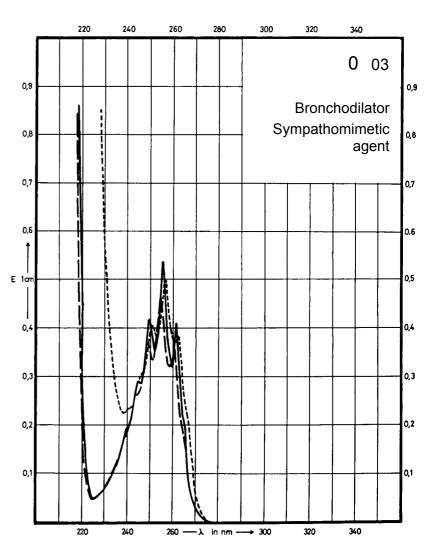

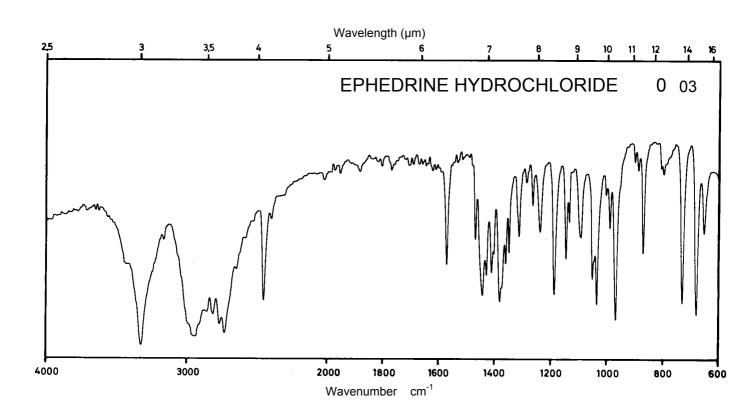





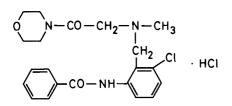


438.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm			
E ^{1%} 1cm	220			
3	9700			

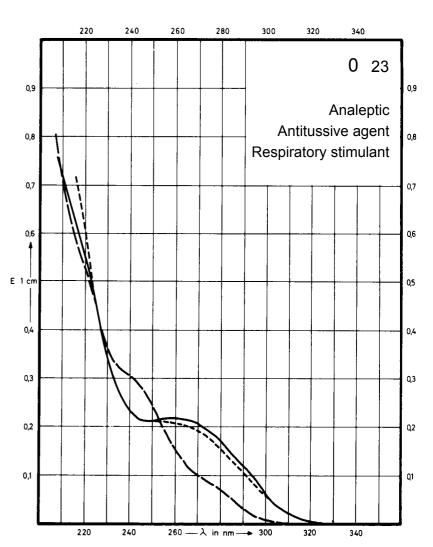

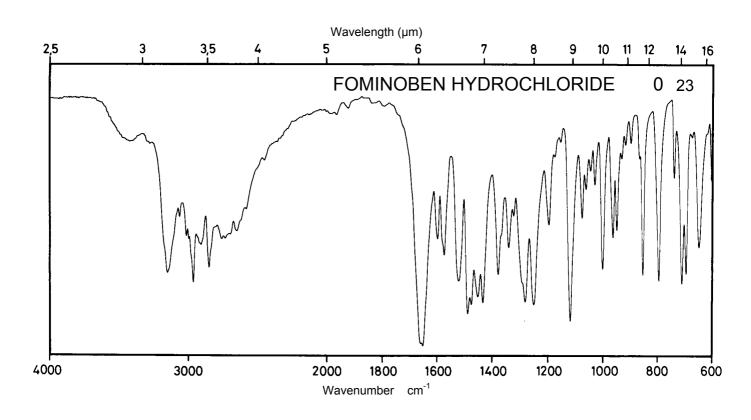




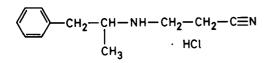



 M_r

224.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm	258 nm		
E ^{1%} 1cm	8.0	8.2		
ε	180	184		

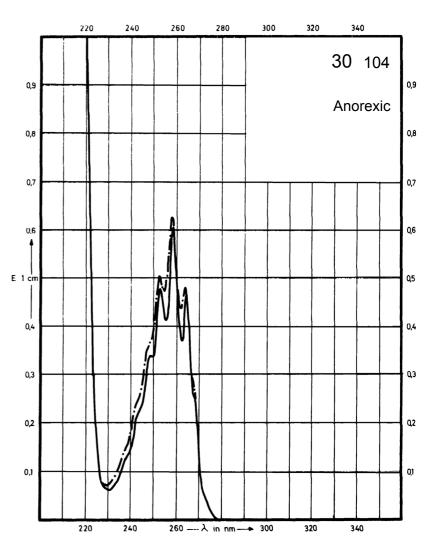


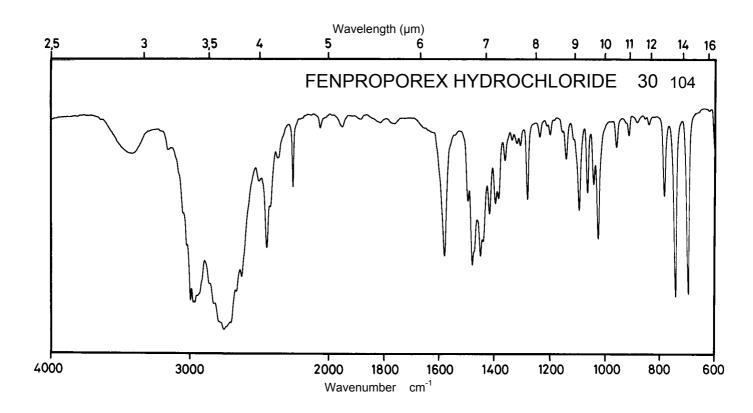


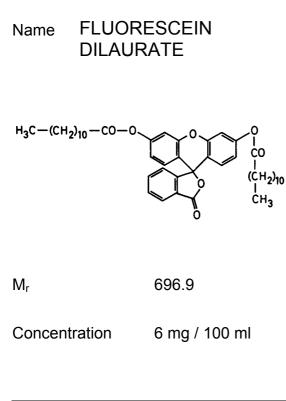


Solvent Symbol	Methanol*	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 282 nm			
E ^{1%} 1cm	74 76			
ε	5160 5300			

^{*} Methanol - Dichlormethane 1 : 1

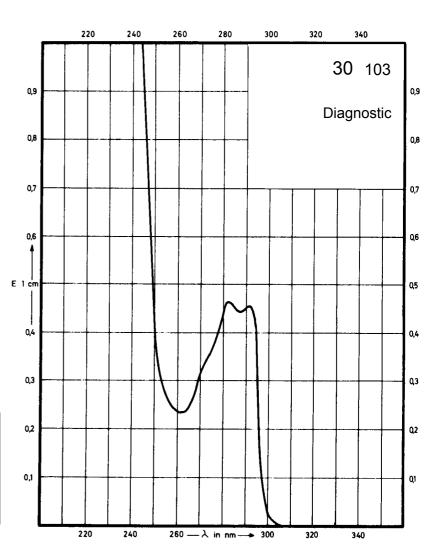

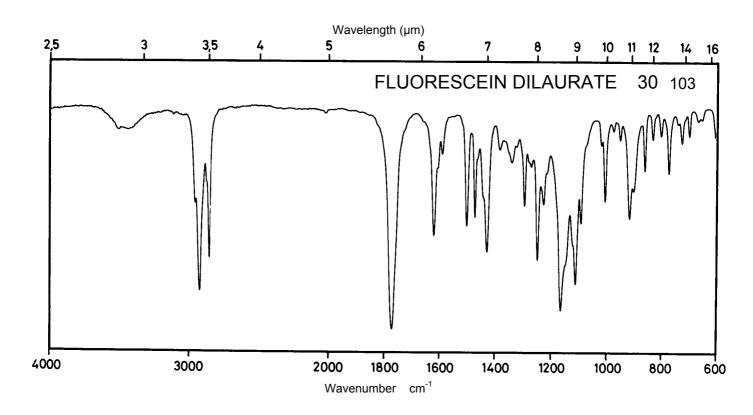




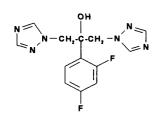


306.3

Concentration 1 mg / 100 ml 26 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm 261 nm			
E ^{1%} 1cm	21.0 23.4			
ε	650 720			

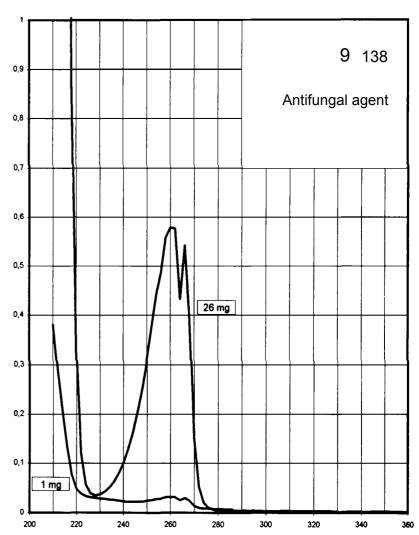

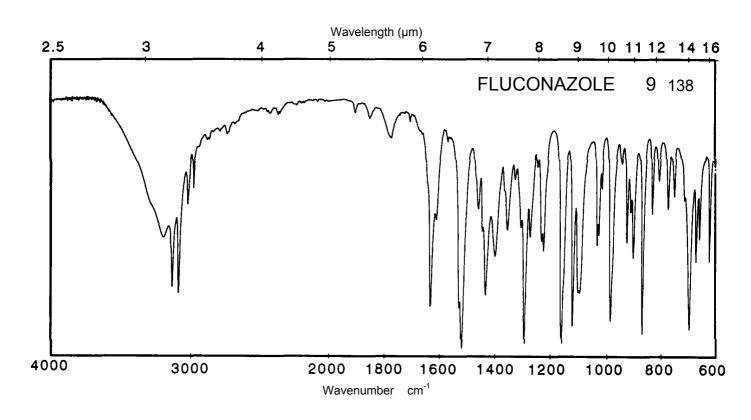



Name FENTICONAZOLE NITRATE

 M_{r}


518.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	253 nm		252 nm	
E ^{1%} 1cm	300		233	
3	15600		12100	

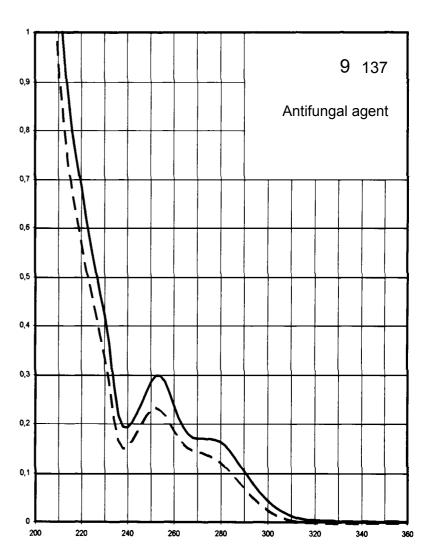

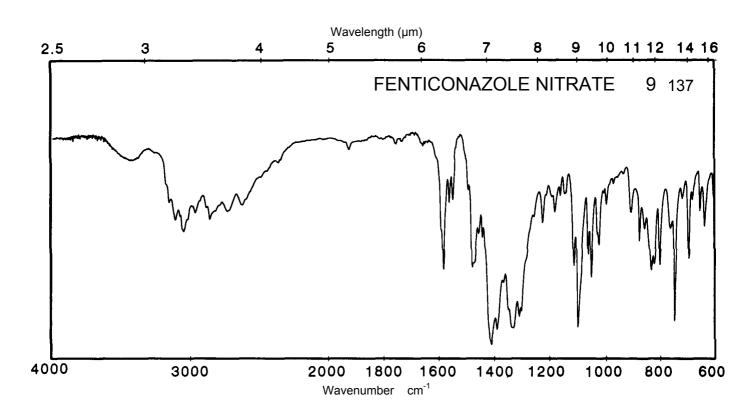



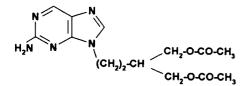



321.3

Concentration 1 n


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	305 nm 243 nm 221 nm	314 nm 225 nm	305 nm 242 nm
E ^{1%} 1cm		221 134 845	131 1234	226 139
ε		7100 4300 27100	4200 39600	7250 4450

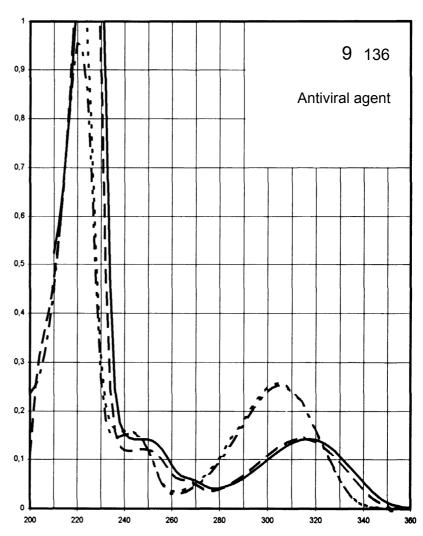

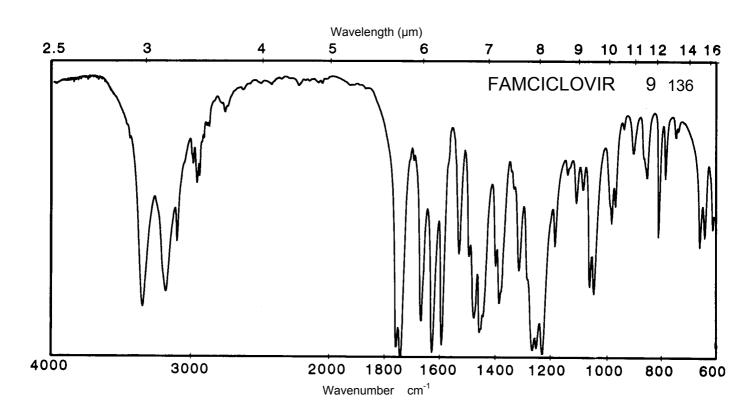





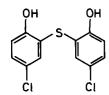


287.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	303 nm		298 nm	323 nm
E ^{1%} 1cm	286		264	352
ε	8200		7600	10100

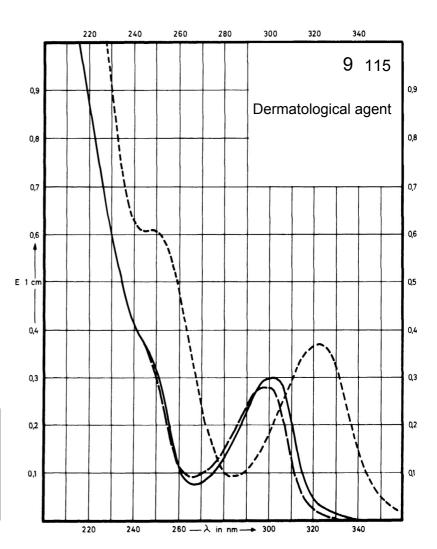

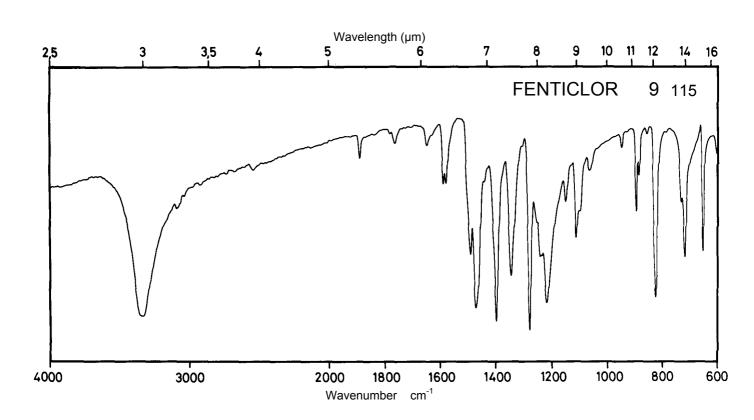




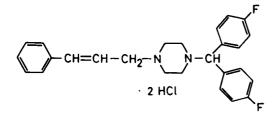


 M_{r} 477.4

Concentration 2 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	293 nm 283 nm 255 nm		254 nm	
E ^{1%} 1cm	23 37 428		421	
ε	1090 1770 20400		20100	

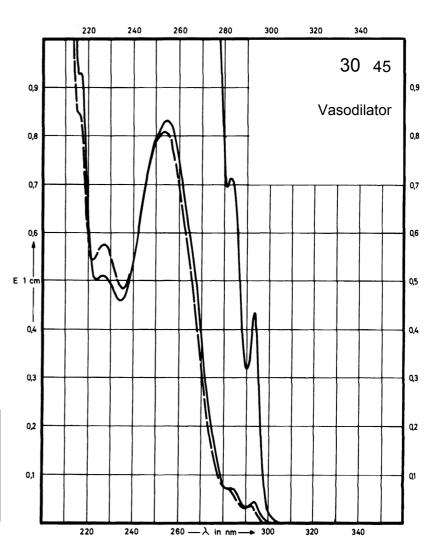

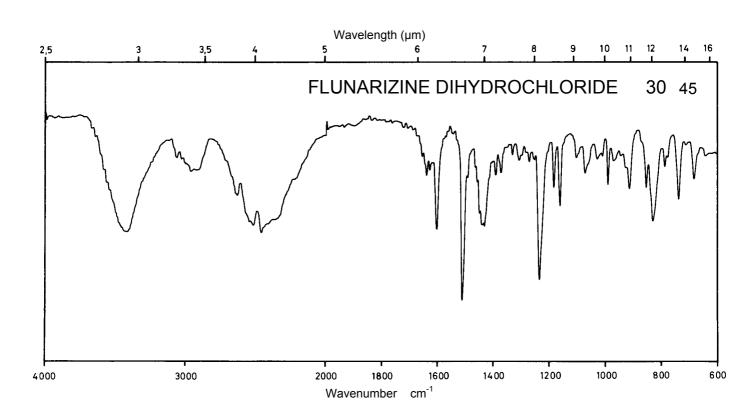


Name FENFLURAMINE HYDROCHLORIDE


$$\begin{array}{c|c} \mathsf{F_3C} & \mathsf{CH_2-CH-CH_3} \\ & \mathsf{I} \\ & \mathsf{NH-CH_2-CH_3} \end{array}$$

 M_{r}


267.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 265 nm	271 nm 265 nm	271 nm 265 nm	
E ^{1%} 1cm	20.0 21.9	20.0 22.5	20.0 22.5	
ε	535 586	535 600	535 600	

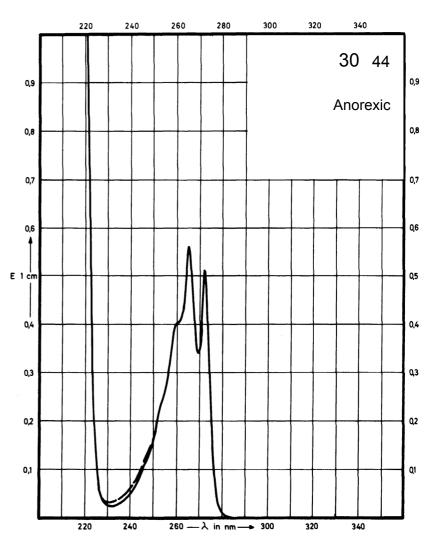

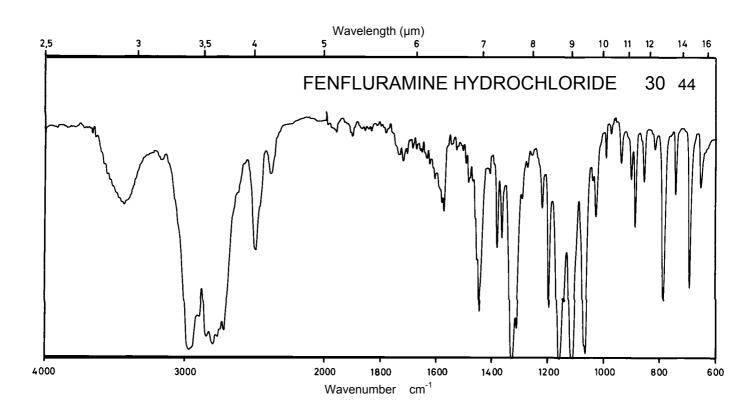





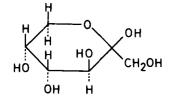

 M_r 180.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

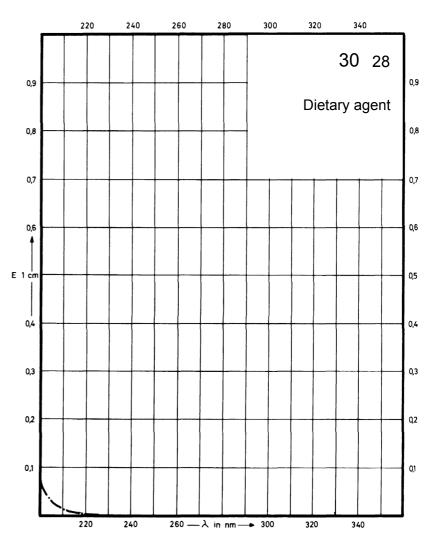

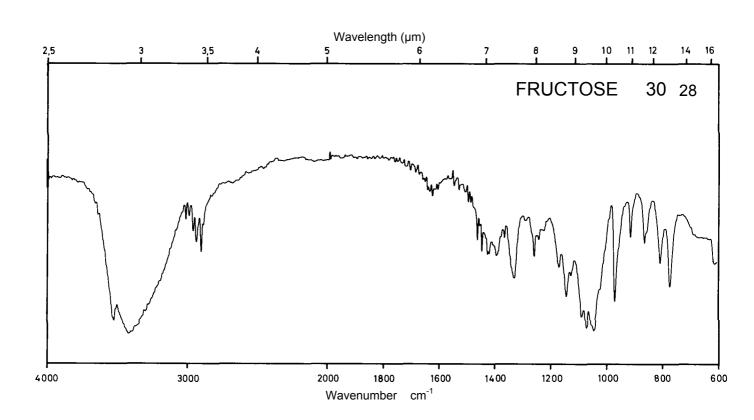

Name FURSULTIAMINE


 M_{r}


398.6

Concentration 2.1


2.1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 234 nm		242 nm	276 nm
E ^{1%} 1cm	141 320		292	141
ε	5610 12760		11630	5640

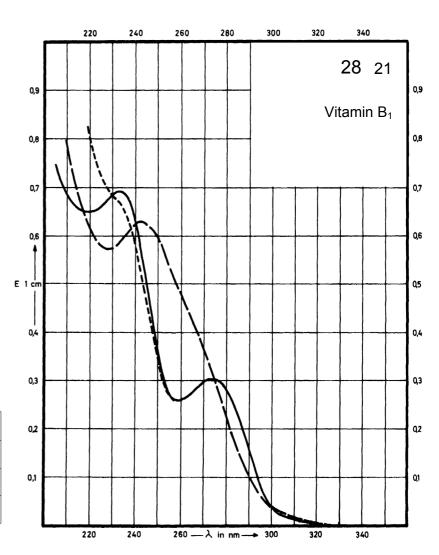

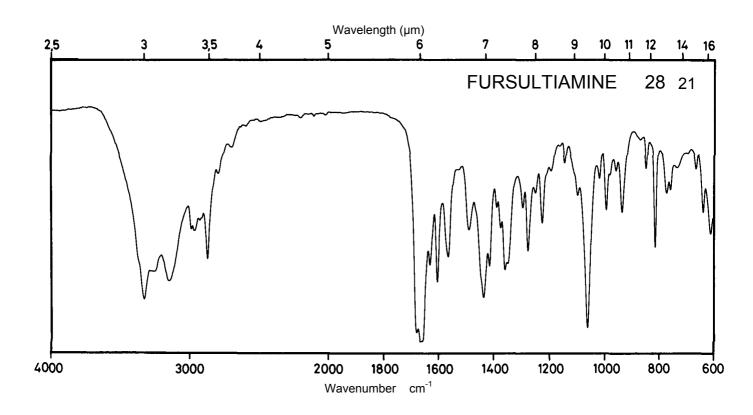


Name FOLIC ACID

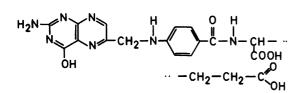


 M_{r}


441.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 260 nm			284 nm 255 nm
E ^{1%} 1cm	657 533			540 551
3	29000 23530			23840 24320

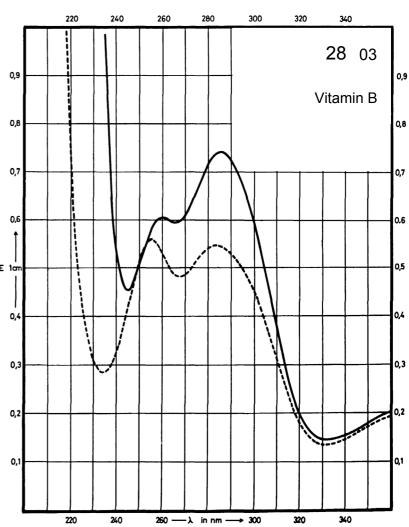

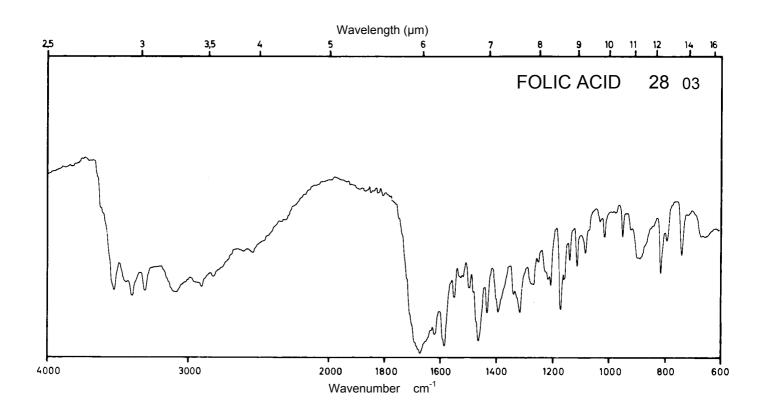





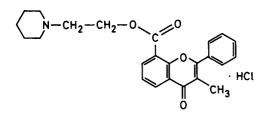


427.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 289 nm 240 nm		319 nm 293 nm 239 nm	Decom- position observed
E ^{1%} 1cm	306 308 419		344 315 388	
ε	13100 13200 17900		14700 13500 16600	

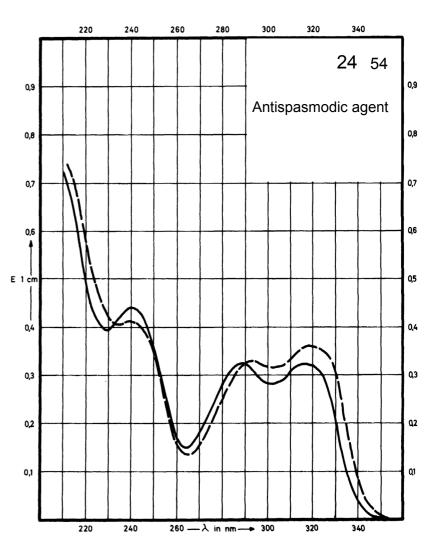

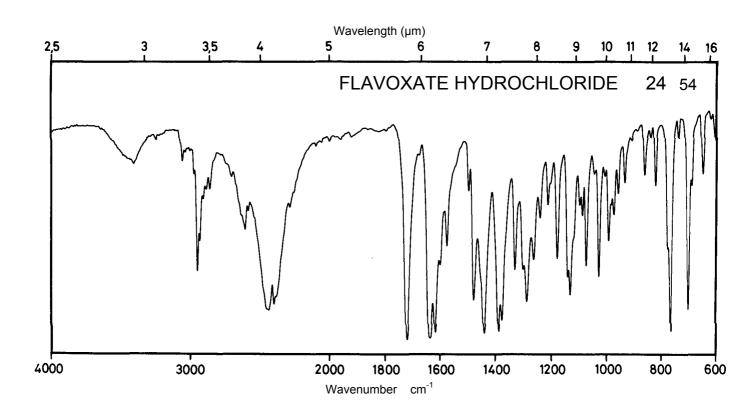



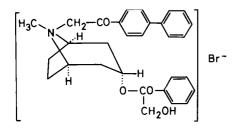



564.5

Concentration


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	297 nm		296 nm	292 nm 255 nm
E ^{1%} 1cm	439		420	265 284
ε	24800		23700	14900 16000

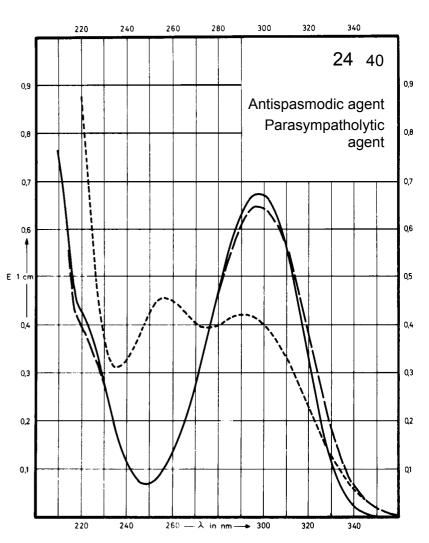

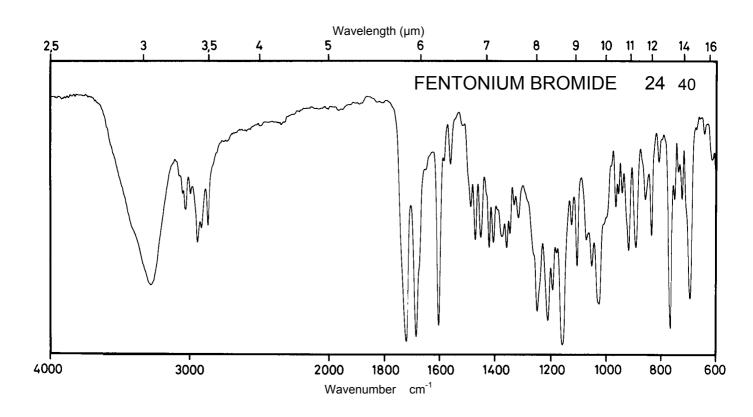


Name FENYRAMIDOL HYDROCHLORIDE

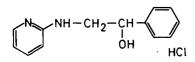


 M_{r}

250.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm 241 nm		311 nm 238 nm	301 nm 241 nm
E ^{1%} 1cm	249 600		265 530	156 576
ε	6240 15030		6630 13290	3920 14440

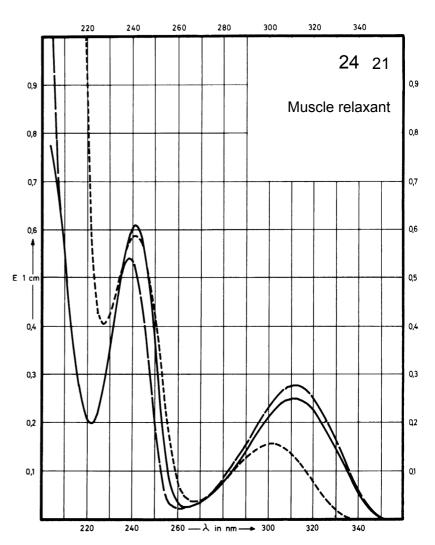

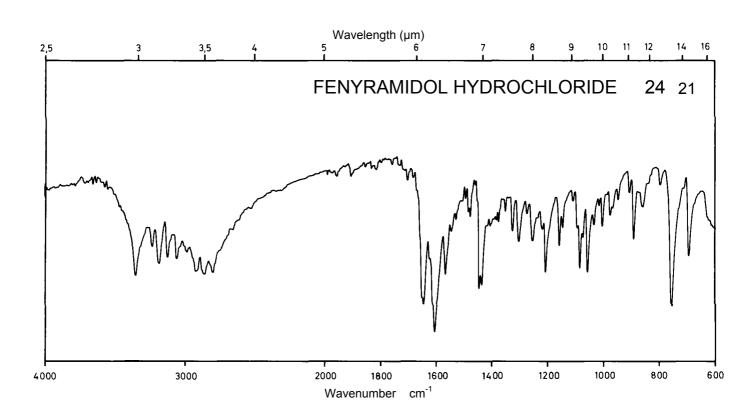



Name FENPIPRANE HYDROCHLORIDE

 M_{r}


315.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm	258 nm	258 nm	
E ^{1%} 1cm	14.0	14.9	14.9	
ε	442	470	470	

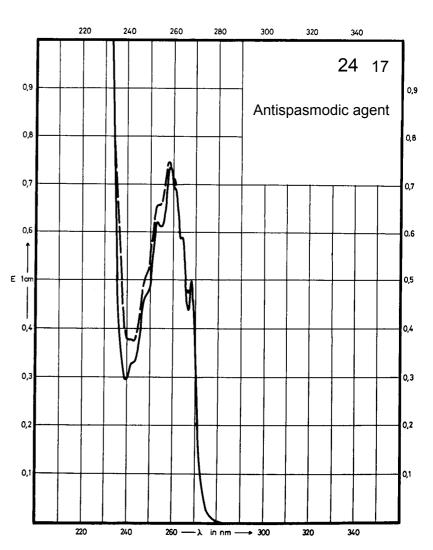

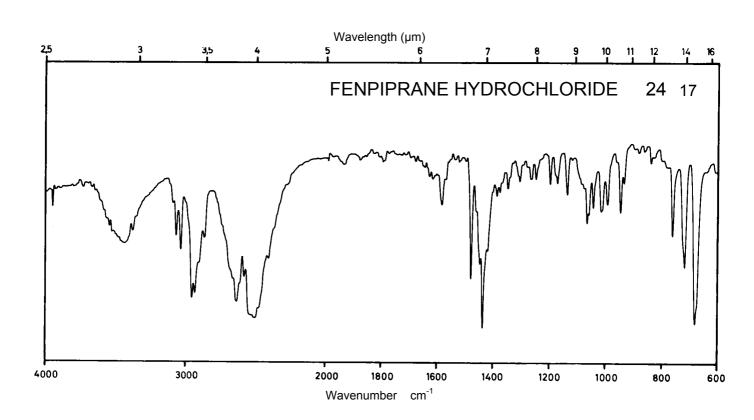

Name FENOTEROL HYDROBROMIDE

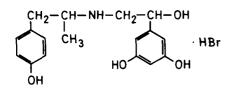



 M_r

384.3

Concentration 2 mg / 100 ml 8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm 224 nm		275 nm 220 nm	295 nm
E ^{1%} 1cm	100 440		86 395	150
ε	3840 16900		3300 15180	5760

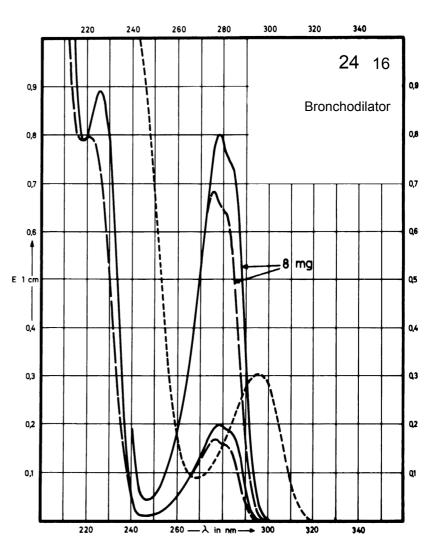

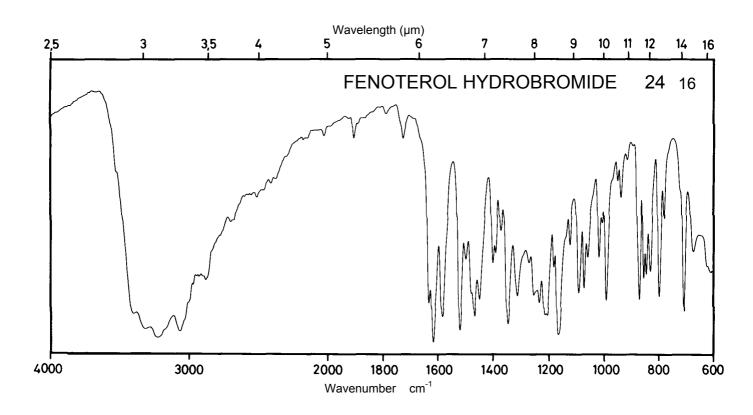





$$\begin{array}{c|c} & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\ & & & \\$$

616.8

Concentration 3 mg / 100 ml 30 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	319 nm 298 nm 243 nm		319 nm 297 nm 286 nm	319 nm 297 nm 286 nm
E ^{1%} 1cm	12 78 262		18 77 113	18 77 115
ε	740 4830 16140		1110 4750 6970	1110 4750 7090

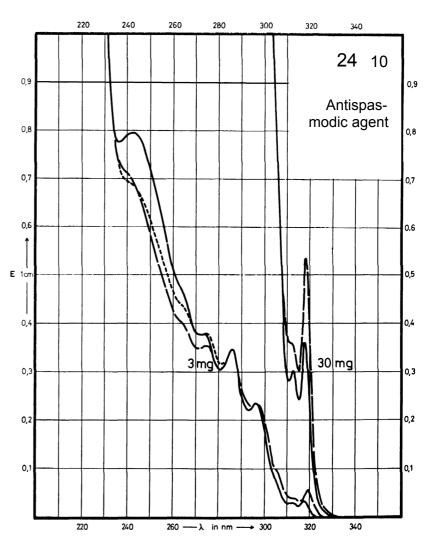

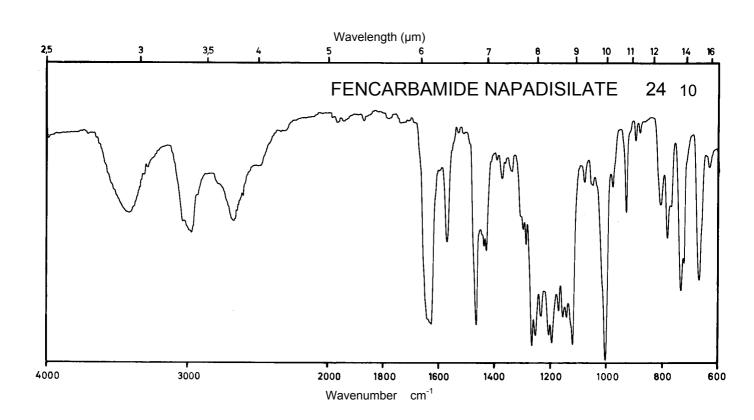


Name FENPIVERINIUM BROMIDE


$$\begin{bmatrix} \downarrow \\ H_2N-CO-C-CH_2-CH_2-N \\ \hline \\ CH_3 \end{bmatrix}^+ Br^-$$

 M_{r}


417.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm		257 nm	258 nm
E ^{1%} 1cm	9.6		10.6	10.4
ε	400		440	430

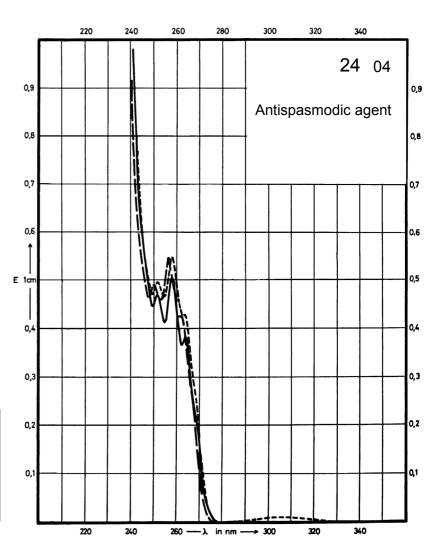

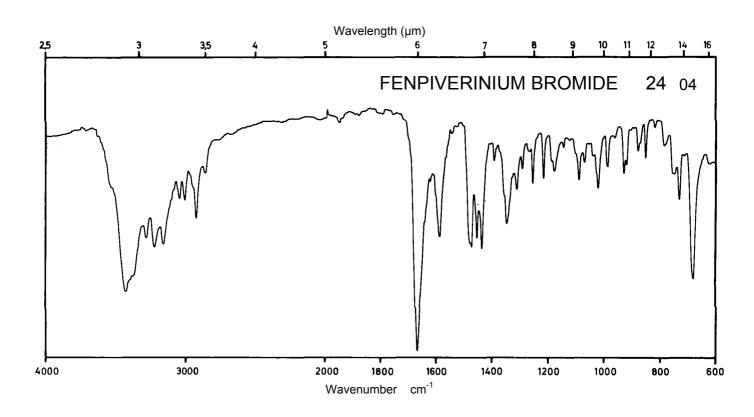





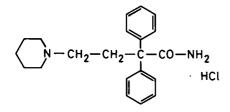


358.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm	258 nm		
E ^{1%} 1cm	11.7	12.5		
ε	420	450		

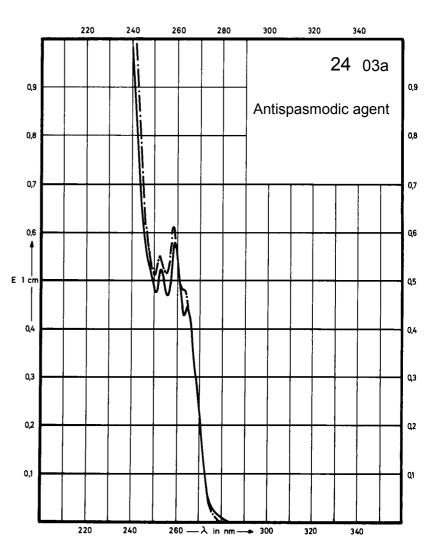

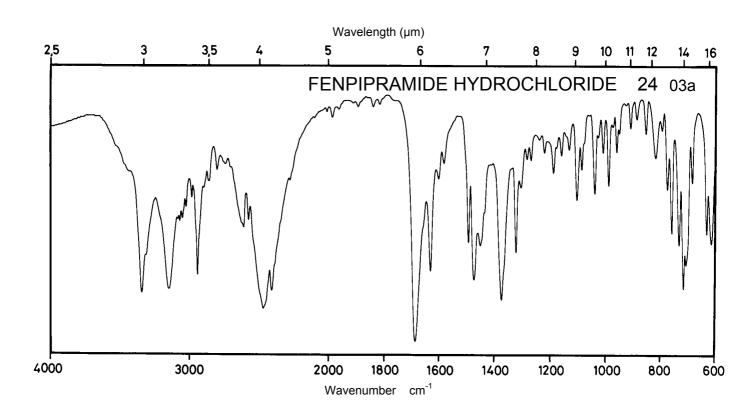




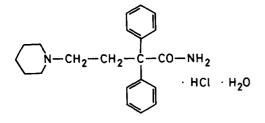



 M_r

376.9


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm		258 nm	
E 1% 1cm	11.3		12.0	
ε	420		450	

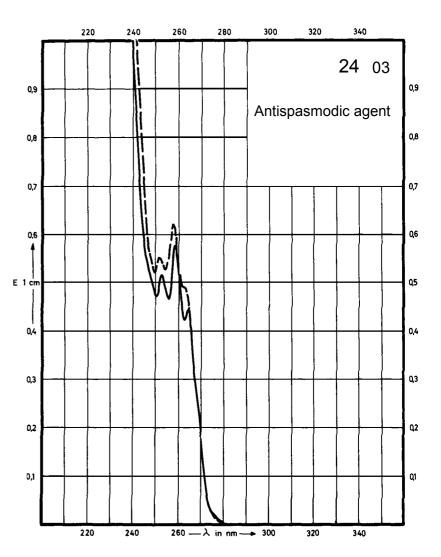


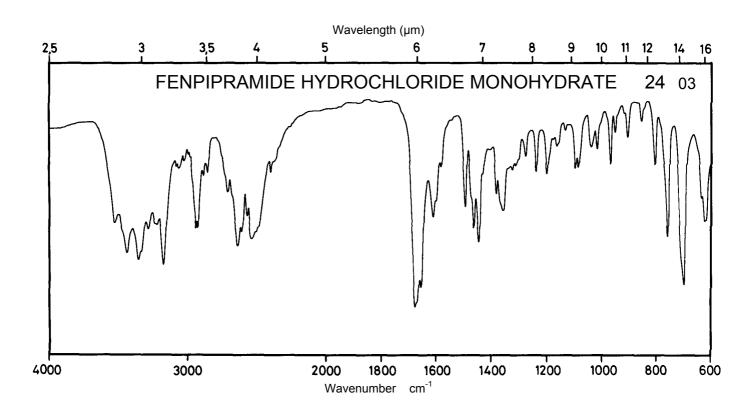

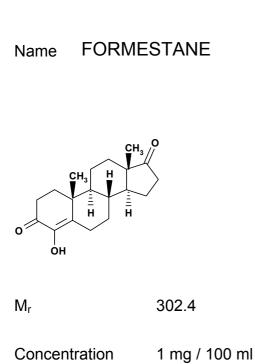


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm	276 nm	276 nm	Decom- position observed
E ^{1%} _{1cm}	408	410	410	
ε	12340	12400	12400	

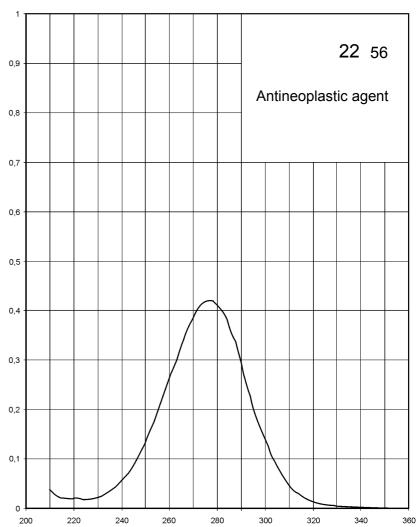

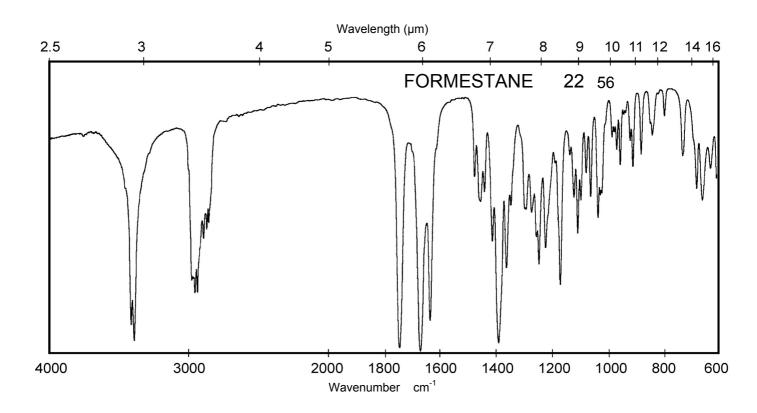





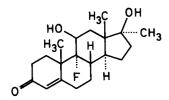


336.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm	243 nm		
E ^{1%} 1cm	499	503		
3	16800	16940		

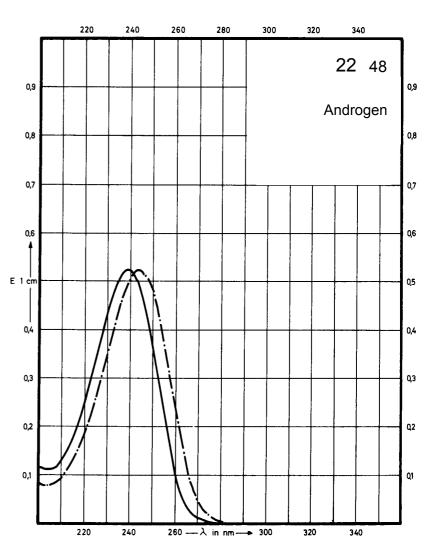

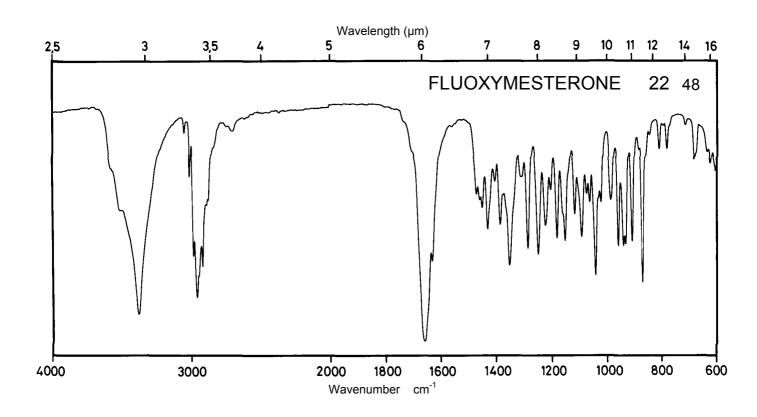




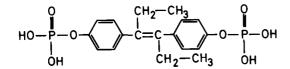


 M_r


428.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	236 nm		240 nm	240 nm
E ^{1%} 1cm	345		328	383
ε	14800		14000	16400

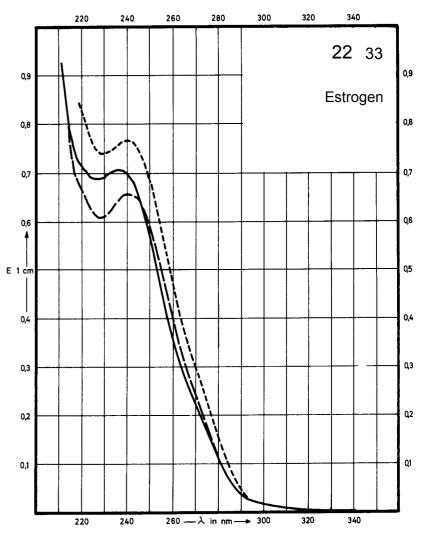

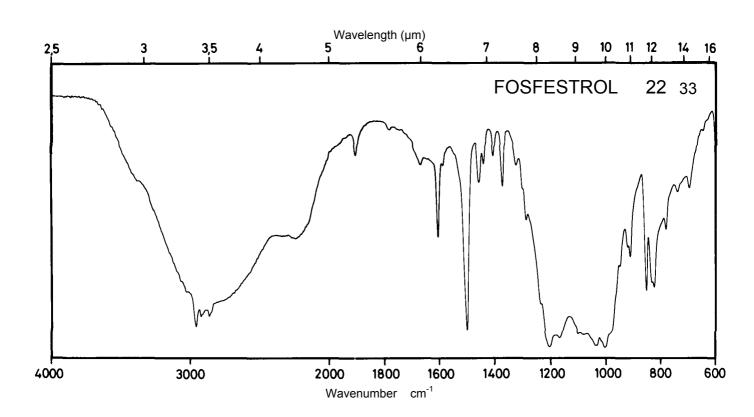





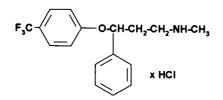

345.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	264 nm 226 nm	263 nm 226 nm	266 nm 225 nm
E ^{1%} 1cm		27.0 365	28.0 367	23.4 354
ε		930 12600	970 12700	810 12200

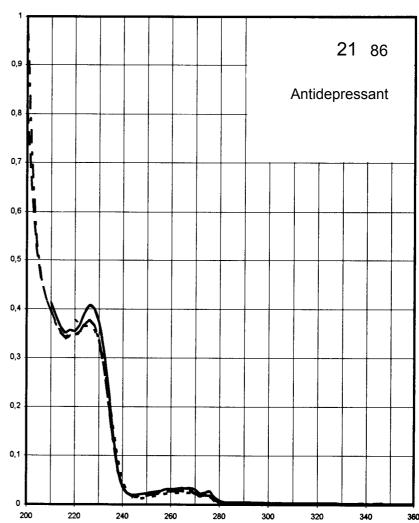

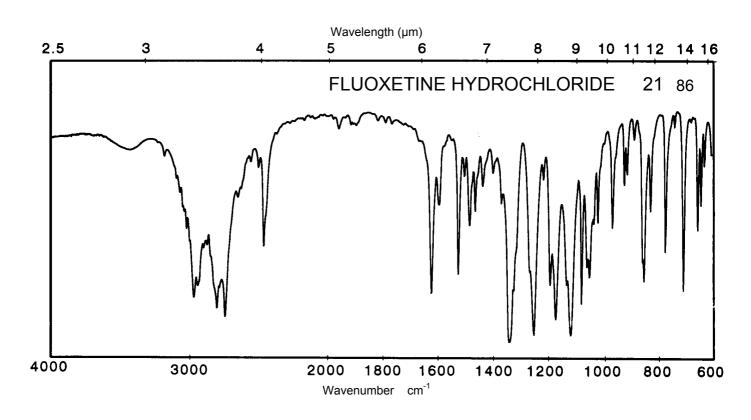



Name FLUVOXAMINE MALEATE

 M_r


434.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	252 nm		246 nm	246 nm
E ^{1%} 1cm	290		258	255
3	12600		11200	11080

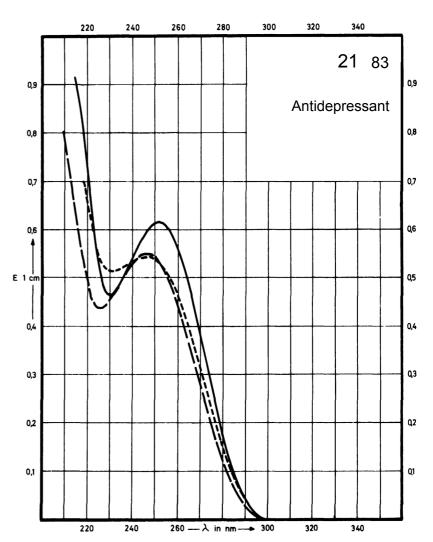

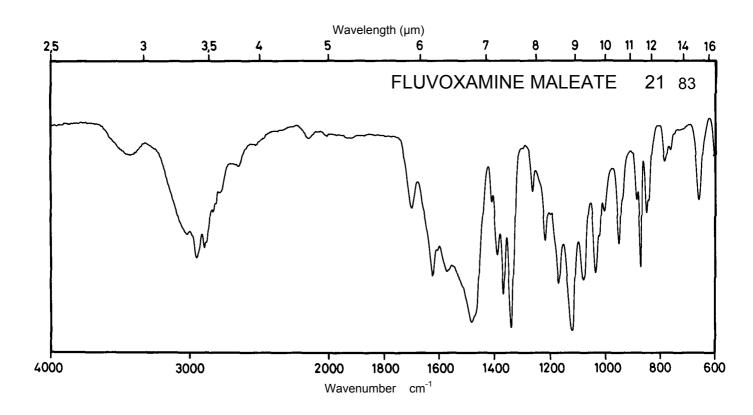



Name FLUPENTIXOL DIHYDROCHLORIDE

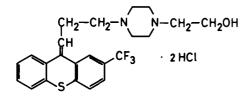



 M_r 507.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm 230 nm	266 nm 230 nm	266 nm 230 nm	266 nm
E ^{1%} 1cm	243 603	233 592	235 610	208
ε	12300 30600	11800 30000	11900 31000	10600

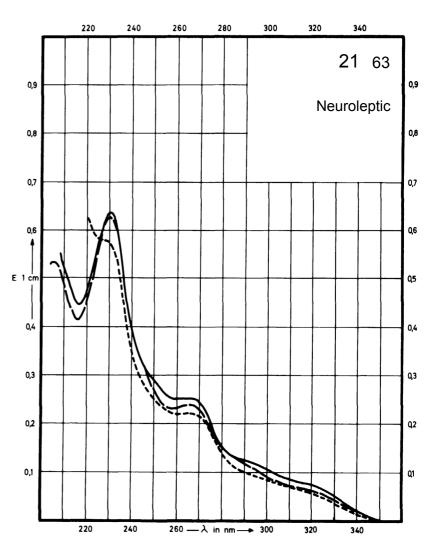


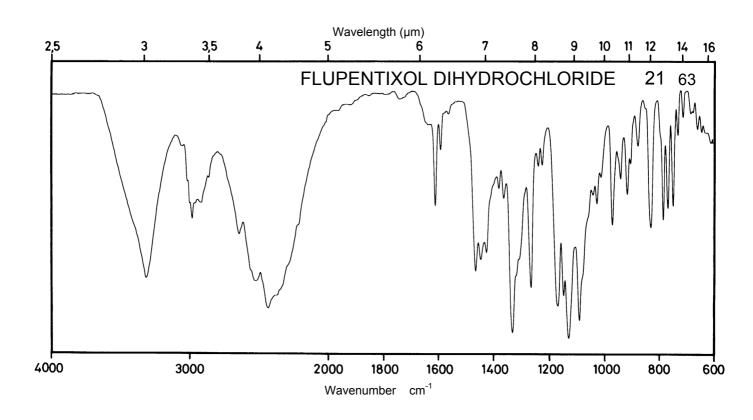

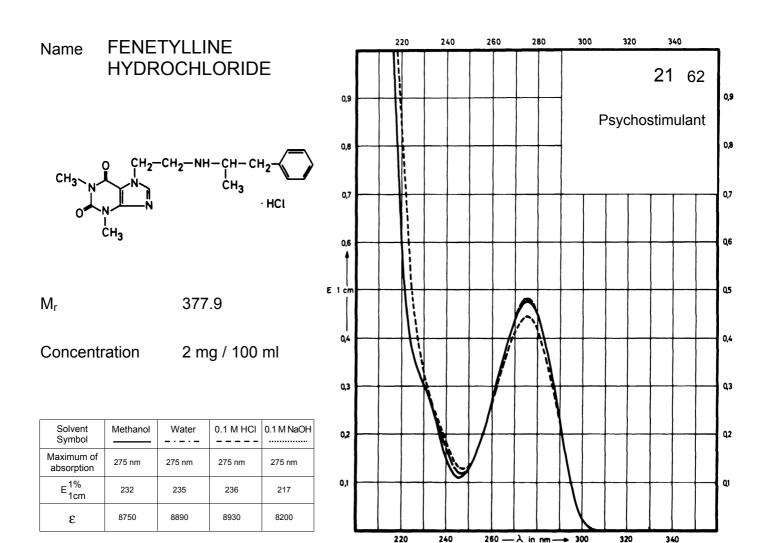

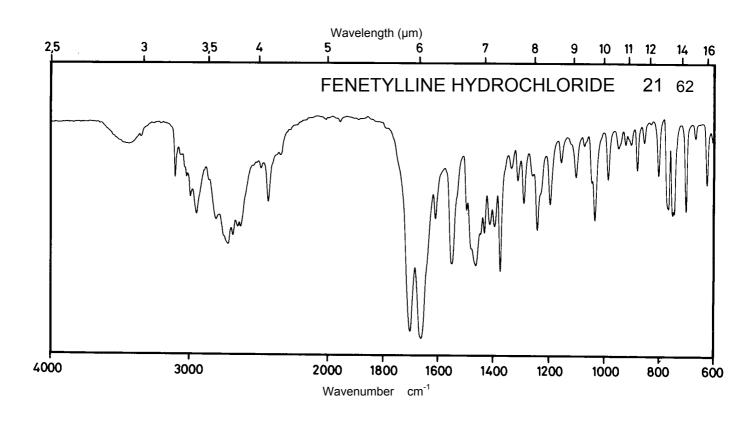



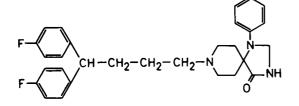



 M_r 475.6

Concentration 2 mg / 100 ml 8 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 N-NaOH*
Maximum of absorption	293 nm 273 nm 253 nm		250 nm	253 nm
E ^{1%} 1cm	39 62 273		163	270
ε	1870 2960 12980		7770	12840

* 1 M NaOH + Methanol (1+9)

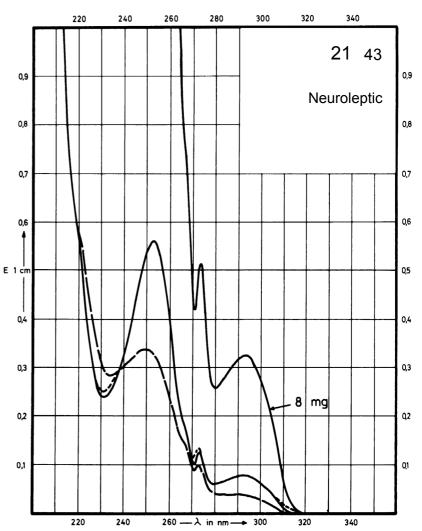

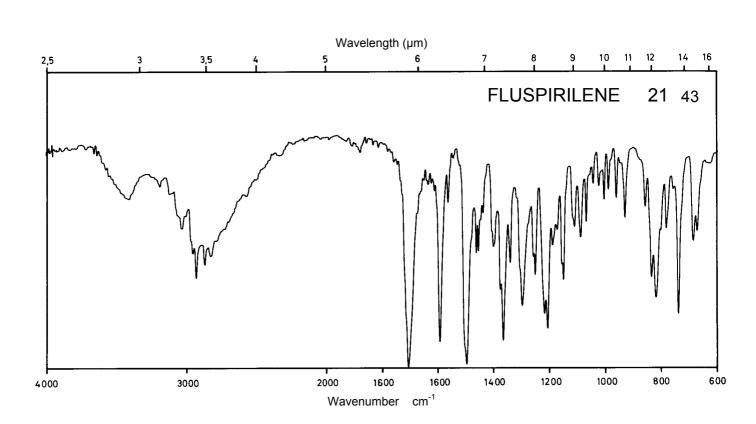

Name FLUPHENAZINE DIHYDROCHLORIDE

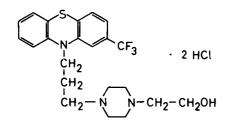



 M_{r}

510.5


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm 260 nm		305 nm 257 nm	309 nm 257 nm
E 1% 1cm	72 651		68 583	69 586
ε	3680 33230		3470 29760	3500 29920

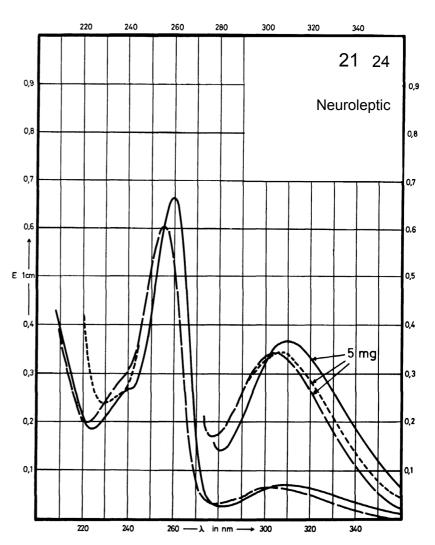

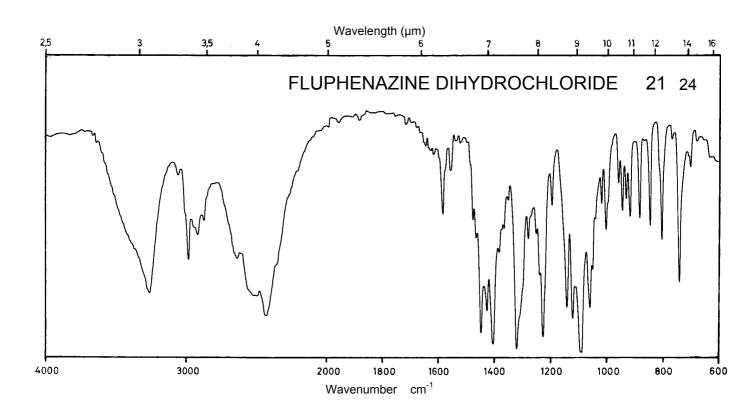





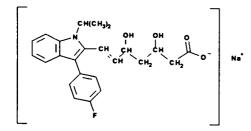



433.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	306 nm 242 nm	Decom- position observed	Decom- position observed	303 nm
E ^{1%} 1cm	294 494			268
ε	12750 21400			11600

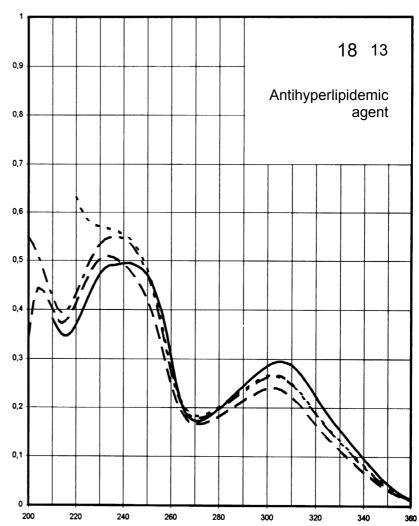


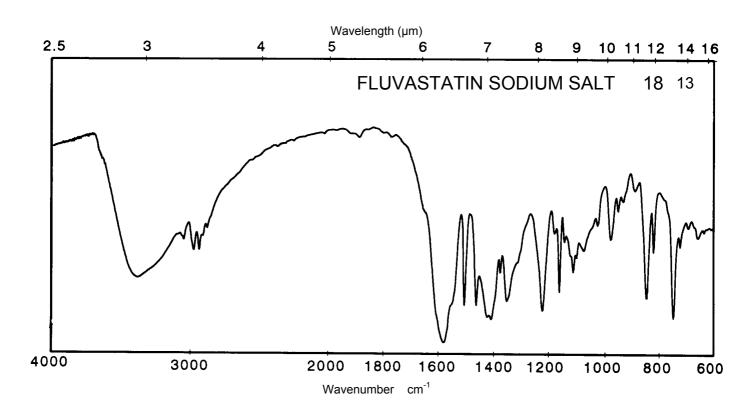

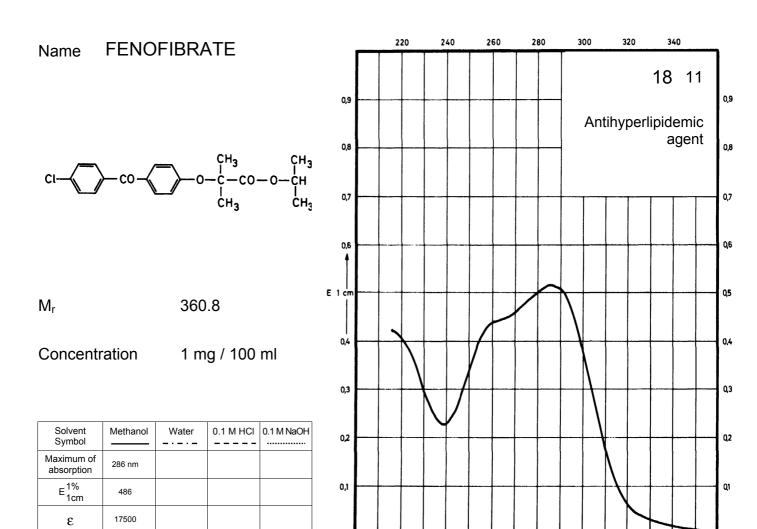

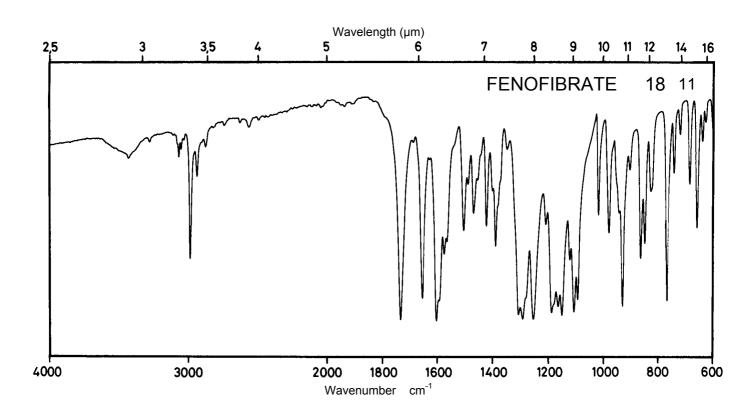



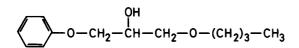



224.3

Concentration


9.4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm 272 nm		276 nm 270 nm	276 nm 270 nm
E ^{1%} _{1cm}	62 74		50 62	50 62
ε	1400 1700		1110 1400	1110 1400

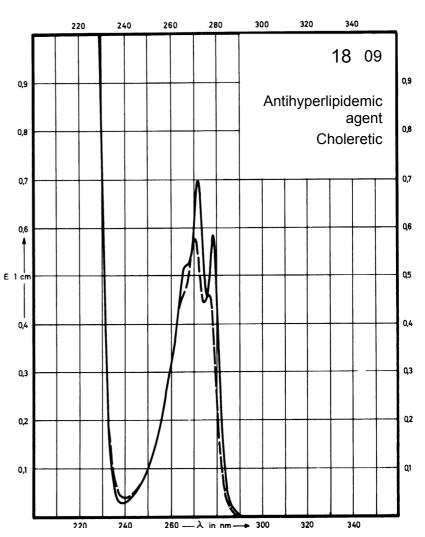

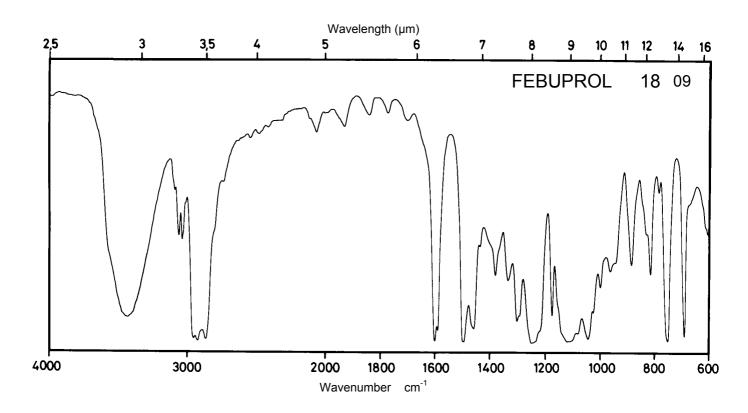





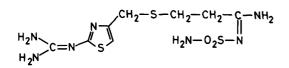


337.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm		265 nm	286 nm
E ^{1%} 1cm	465		309	440
ε	15700		10400	14850

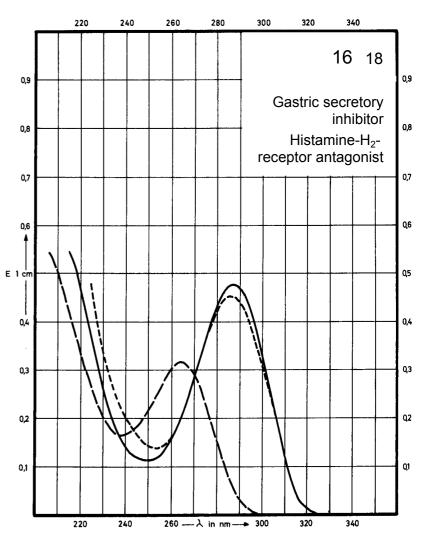

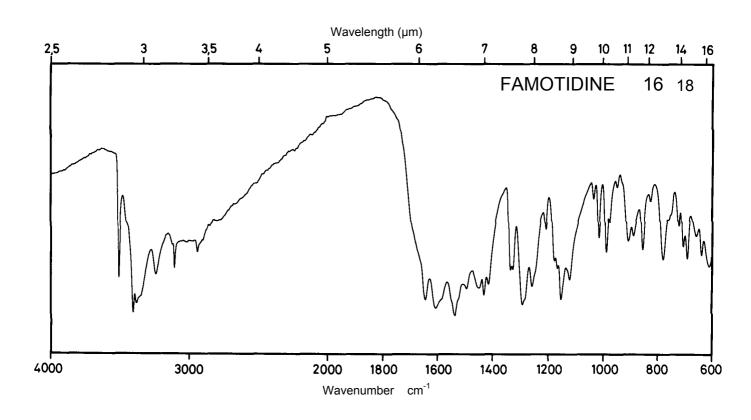

Name FOSAZEPAM

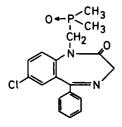



 M_{r}

360.8


Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	313 nm 229 nm	231 nm	357 nm 287 nm 240 nm	Decom- position observed
E ^{1%} 1cm	58 810	855	108 378 708	
ε	2090 29220	30850	3900 13640 25540	

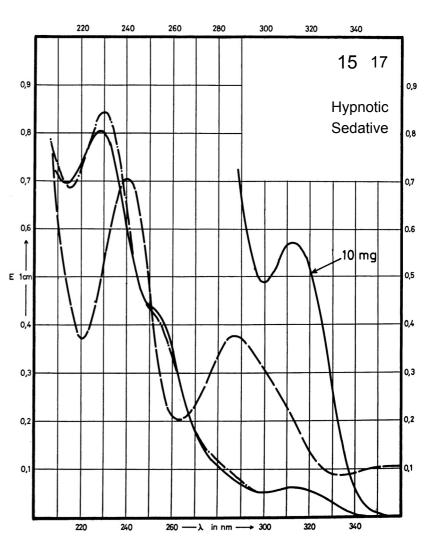

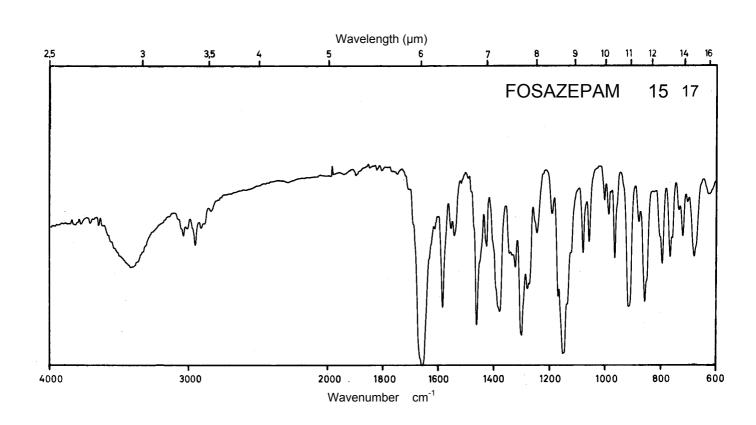


Name FLUNITRAZEPAM

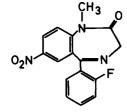


 M_{r}


313.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	308 nm 252 nm		Decom- position observed	Decom- position observed
E ^{1%} 1cm	332 523			
ε	10400 16390			

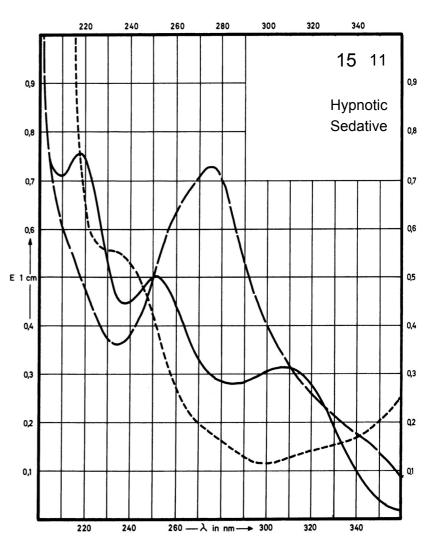

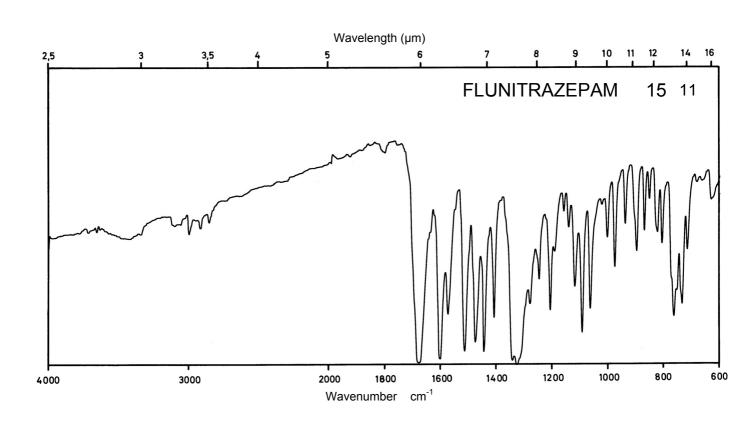


Name FLURAZEPAM DIHYDROCHLORIDE

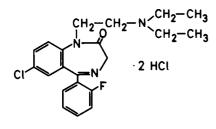


 M_r 460.8


Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	312 nm 228 nm			314 nm 230 nm
E ^{1%} 1cm	43 742			45 721
3	1970 34220			2060 33220

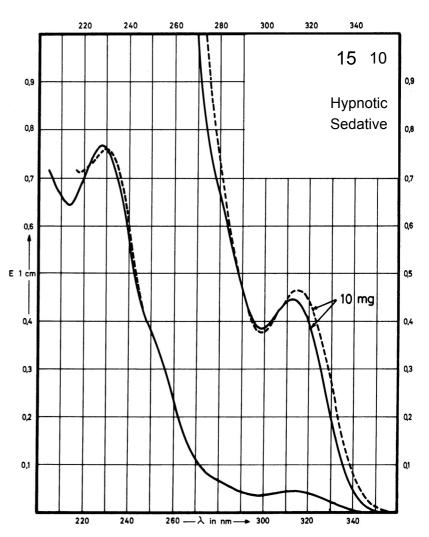

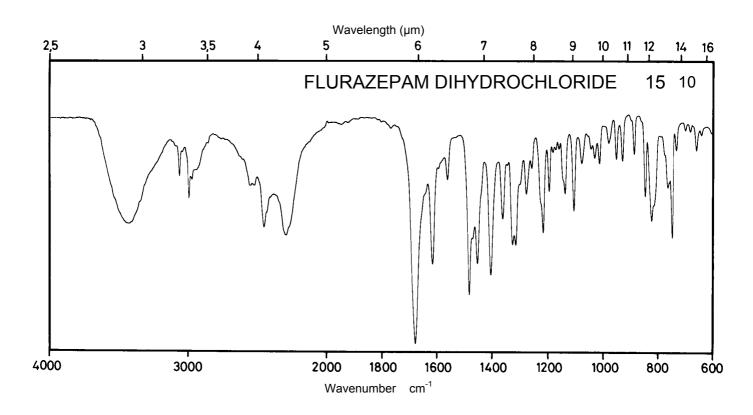





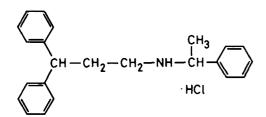


351.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm			
E ^{1%} 1cm	16.8			
ε	590			

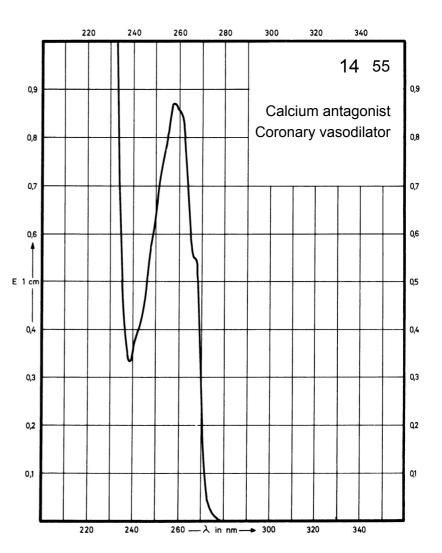

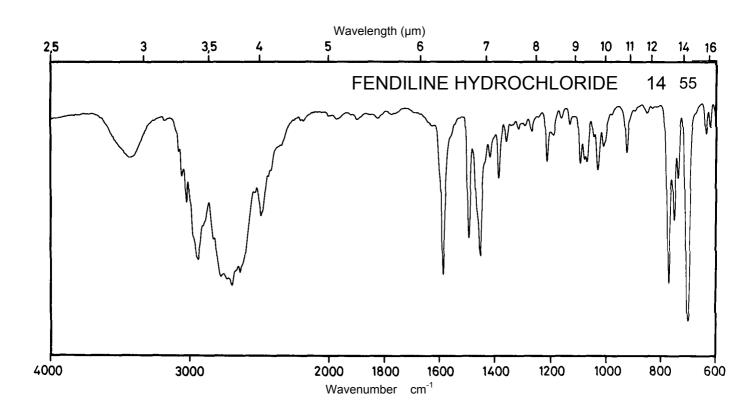




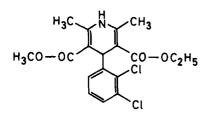


 M_r

384.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	362 nm 237 nm			
E ^{1%} 1cm	186 524			
ε	7150 20100			

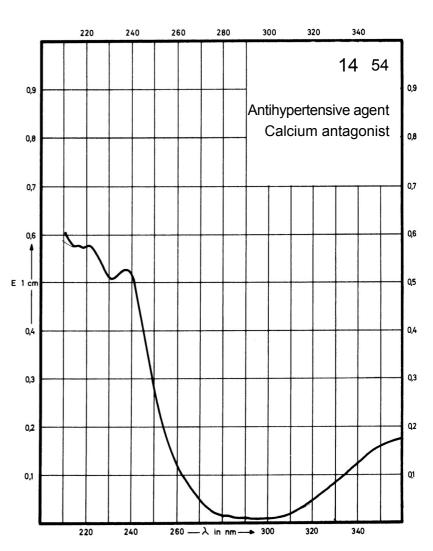

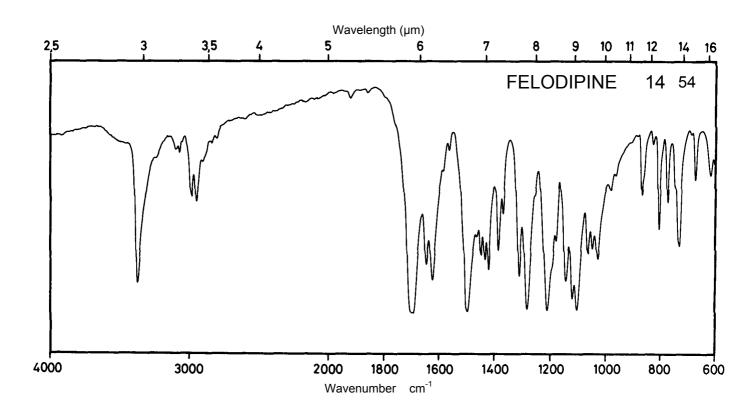



Name FLECAINIDE ACETATE

 M_{r}


474.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	297 nm		296 nm	294 nm
E ^{1%} 1cm	63		60	58
3	3000		2800	2750

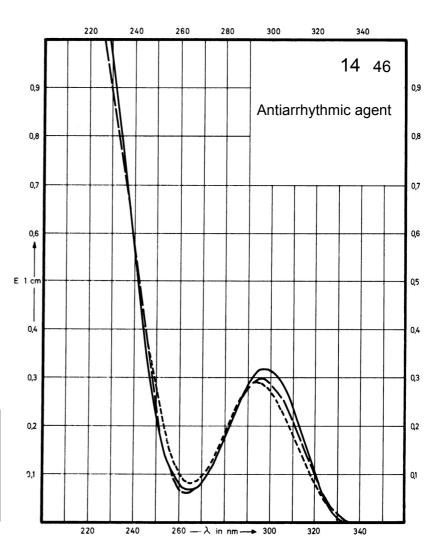

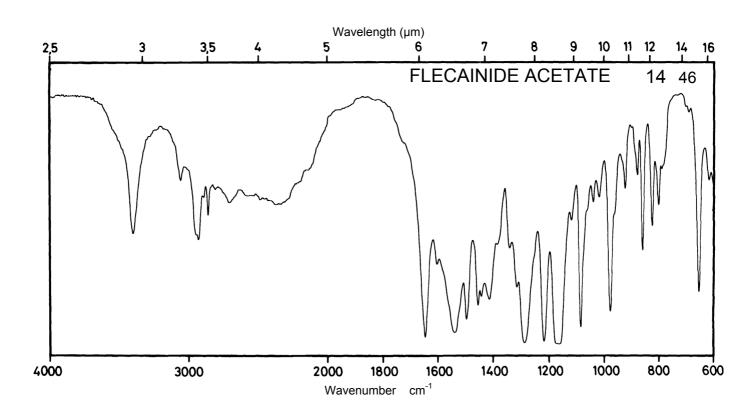



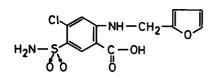



330.7

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	337 nm 273 nm 233 nm		340 nm 274 nm 235 nm	335 nm 270 nm
E ^{1%} _{1cm}	159 671 1390		169 638 1430	131 586
3	5260 22190 45970		5590 21100 47290	4330 19380

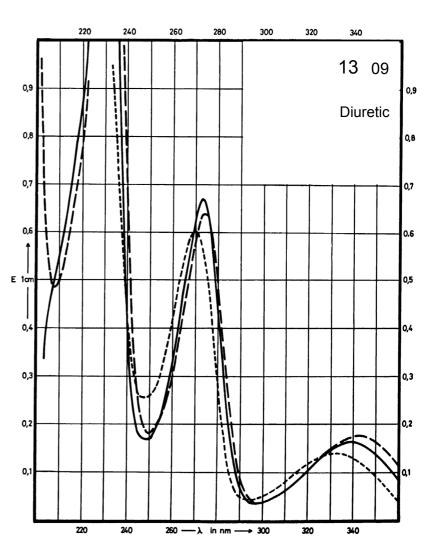

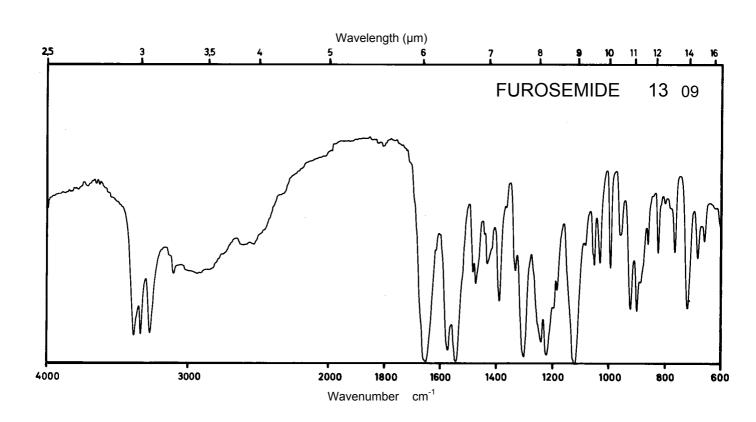





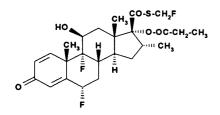


500.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	237 nm			
E ^{1%} 1cm	372			
3	18600			

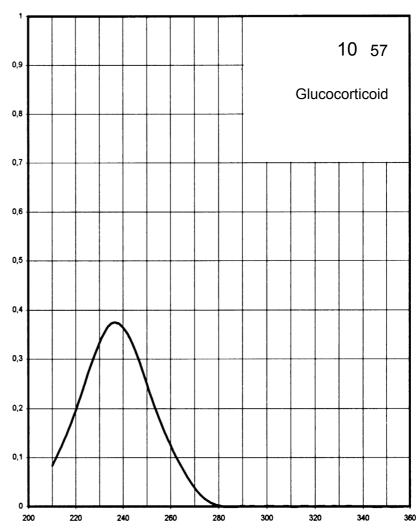

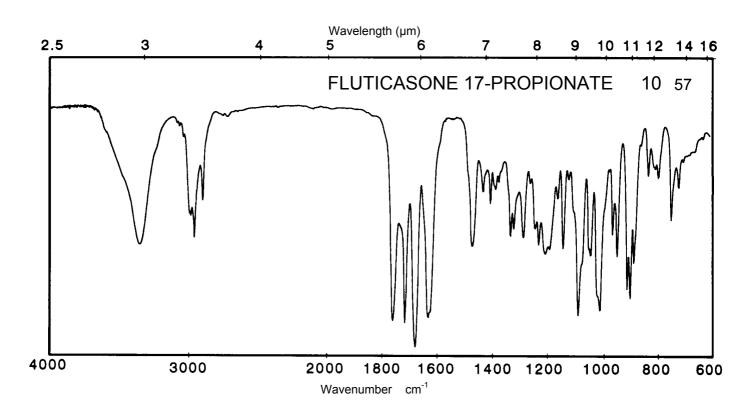





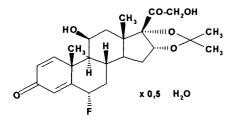


 M_r 443.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	242 nm	245 nm	245 nm	245 nm
E ^{1%} 1cm	370	384	385	376
ε	16400	17000	17000	16700

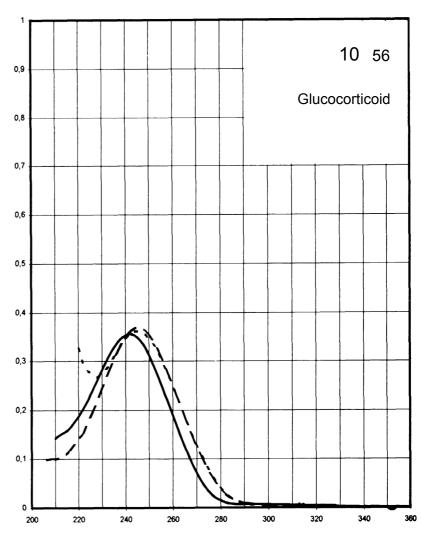

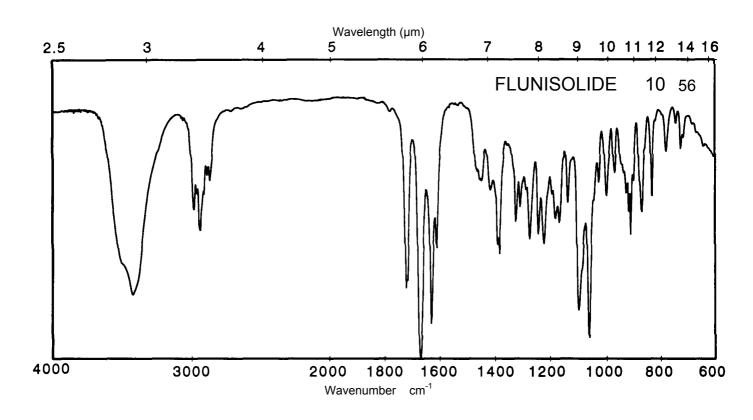

Name FLUPREDNIDENE 21-ACETATE

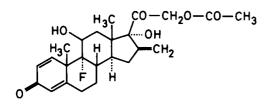



432.5

Concentration


1.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm	241 nm		
E ^{1%} 1cm	361	365		
ε	15600	15800		

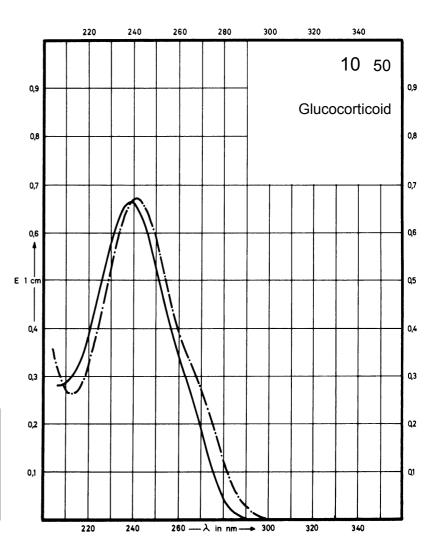

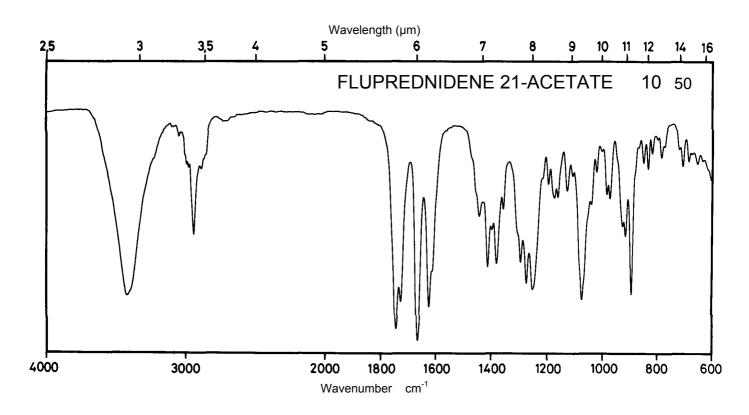





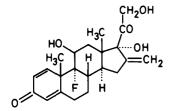


 M_r 390.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm	241 nm		
E ^{1%} 1cm	400	403		
ε	15600	15750		

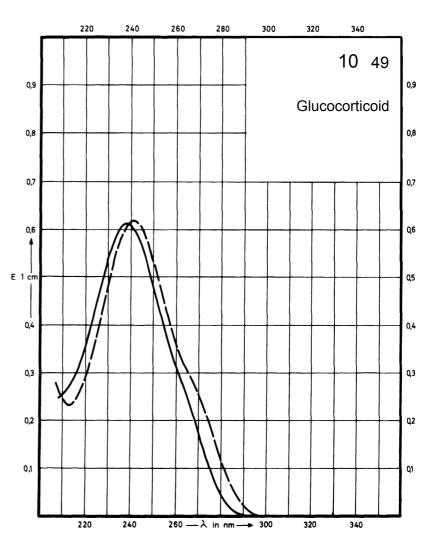

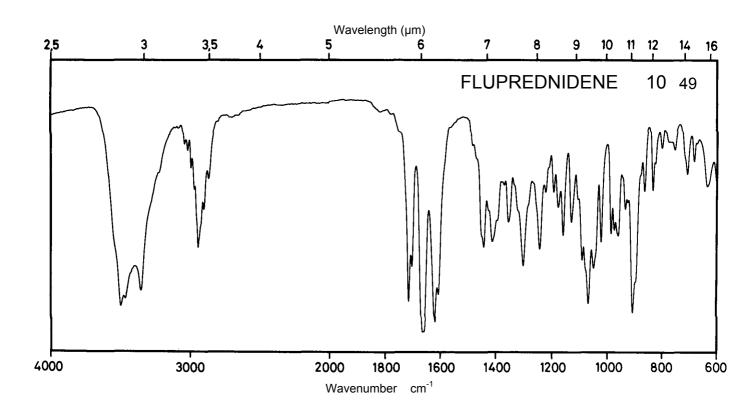





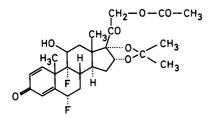

 M_r 494.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	237 nm			
E ^{1%} 1cm	335			
3	16550			

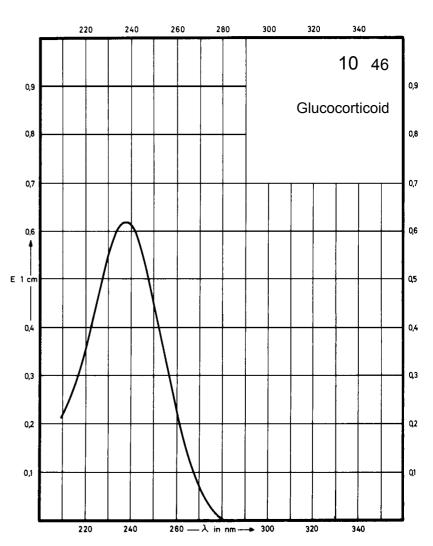

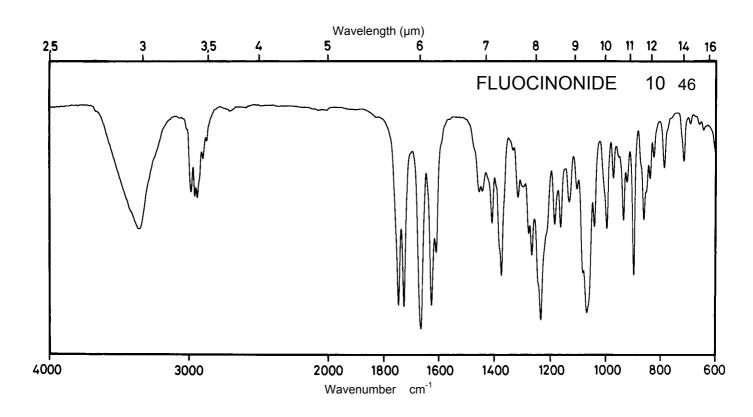



Name FLUOCINOLONE ACETONIDE

 M_{r}


452.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	237 nm	240 nm		
E ^{1%} 1cm	358	363		
3	16200	16400		

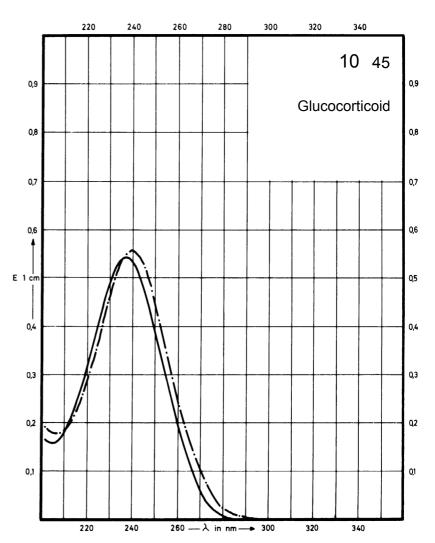

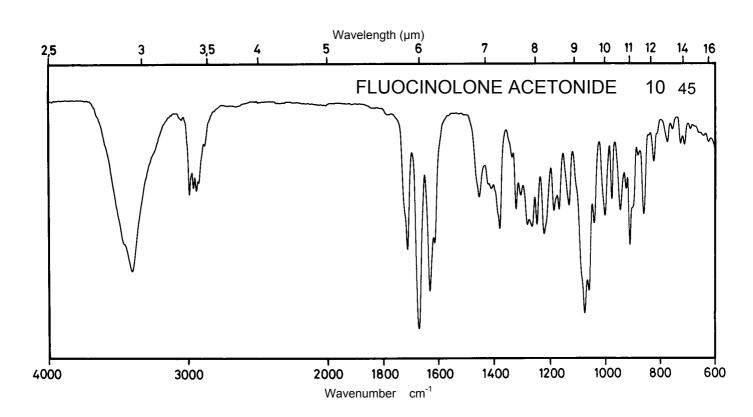





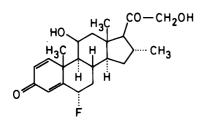

376.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	242 nm	246 nm		
E ^{1%} 1cm	428	440		
3	16120	16570		

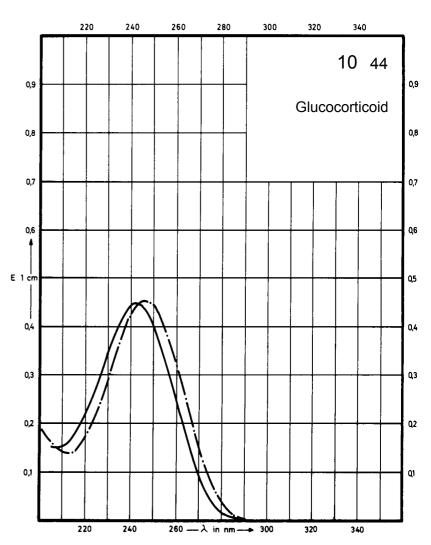

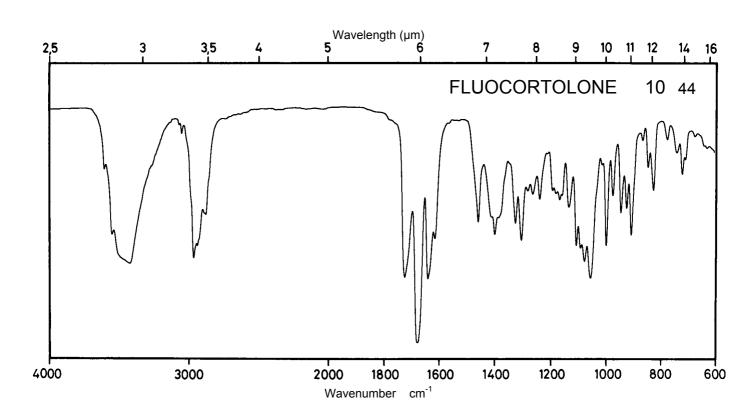



Name FLUOCORTIN BUTYL

 M_{r}


427.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm			
E ^{1%} 1cm	379			
3	16210			

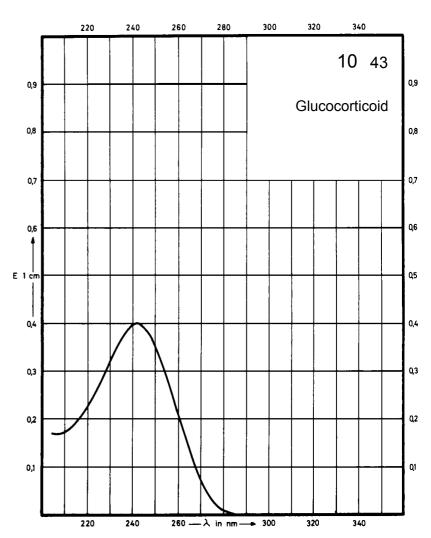

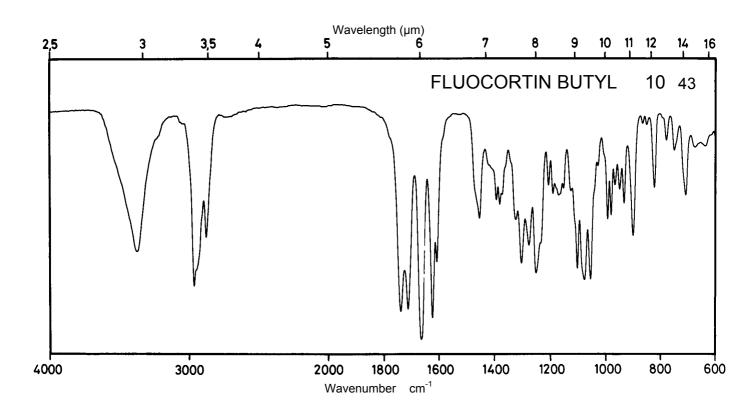



Name FLUDROXYCORTIDE

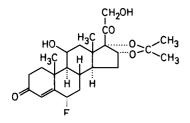


 M_r 436.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	236 nm	241 nm		
E ^{1%} 1cm	350	350		
ε	15300	15300		

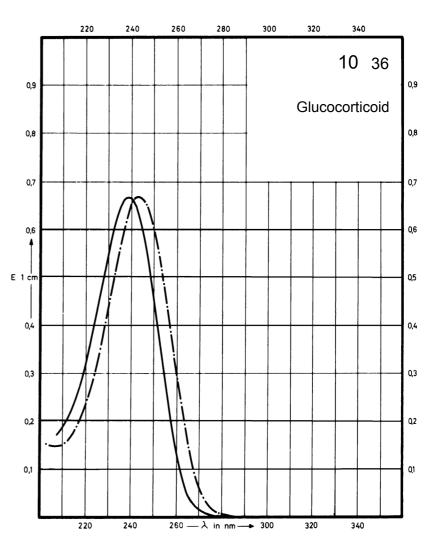

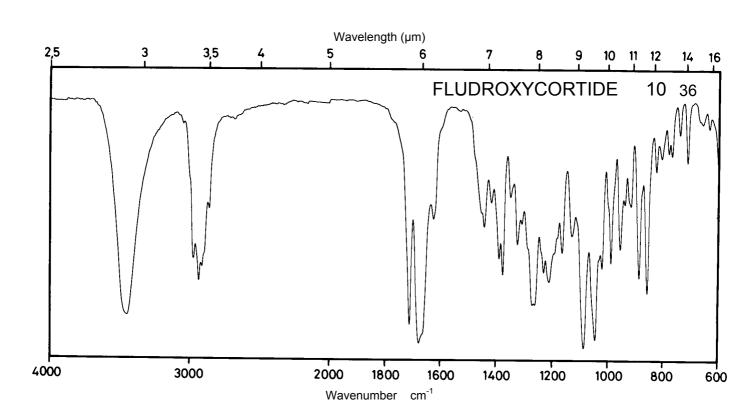



Name FLUDROCORTISONE 21-ACETATE

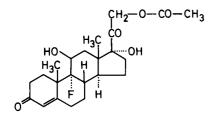

 M_r 422.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm	243 nm		
E ^{1%} 1cm	411	415		
3	17400	17500		

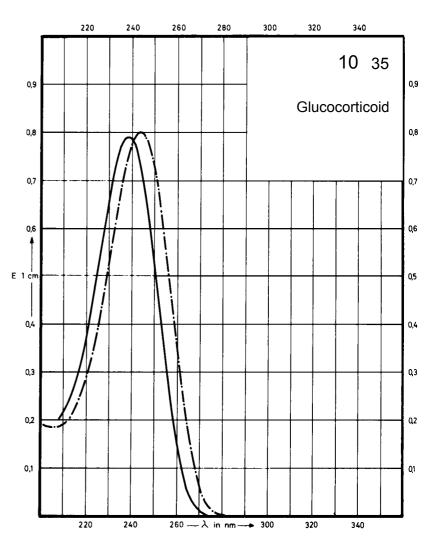


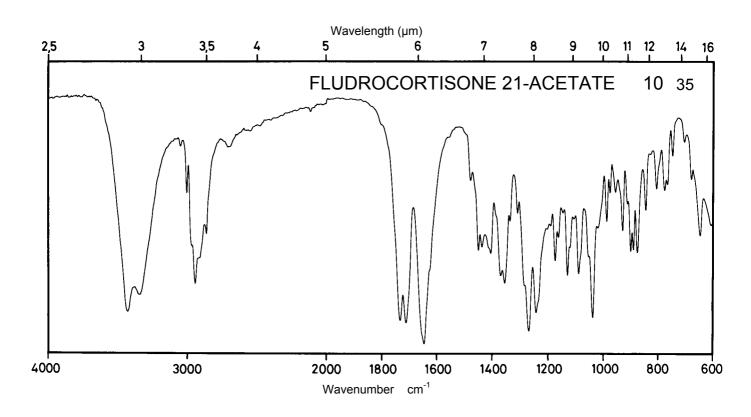


Name FLUDROCORTISONE


 M_r 380.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm	243 nm		
E ^{1%} 1cm	456	456		
ε	17400	17400		

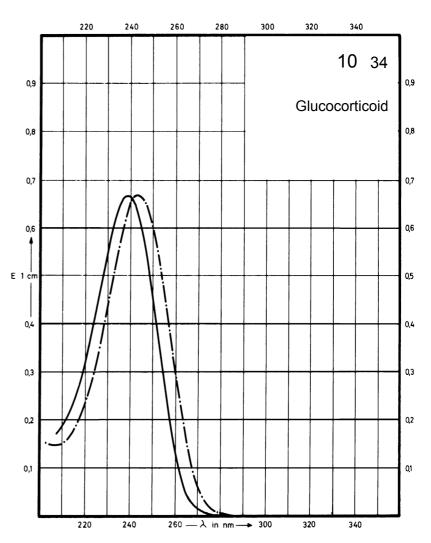

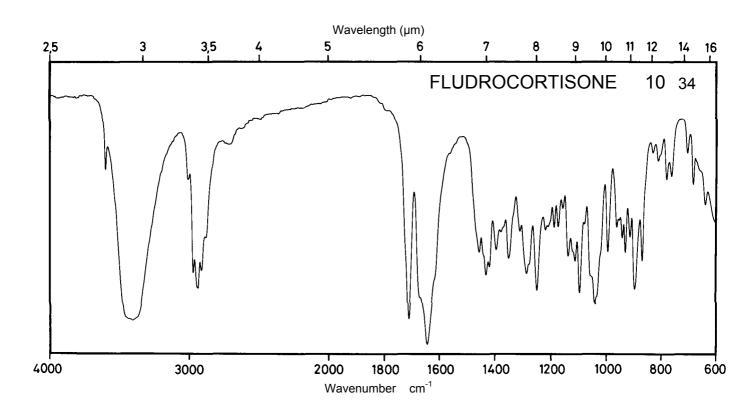



Name FLUMETASONE 21-PIVALATE

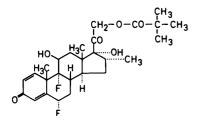


494.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm	242 nm		
E ^{1%} 1cm	328	318		
ε	16230	15720		

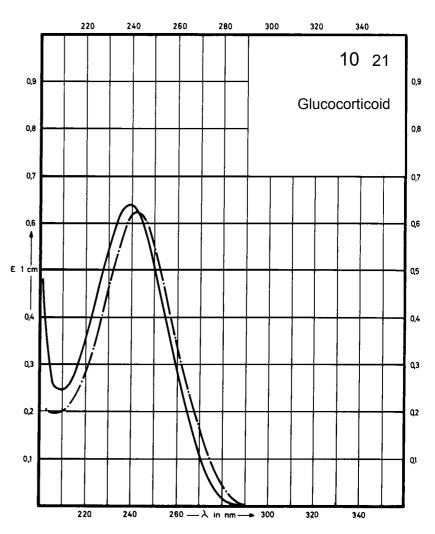

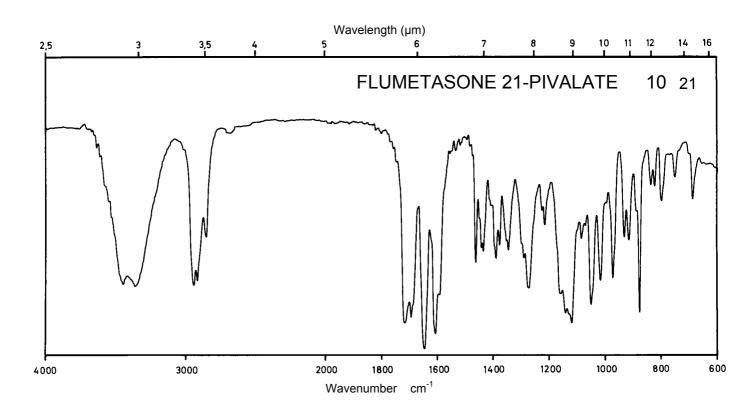





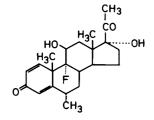


 M_r 376.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm			
E ^{1%} 1cm	410			
3	15440			

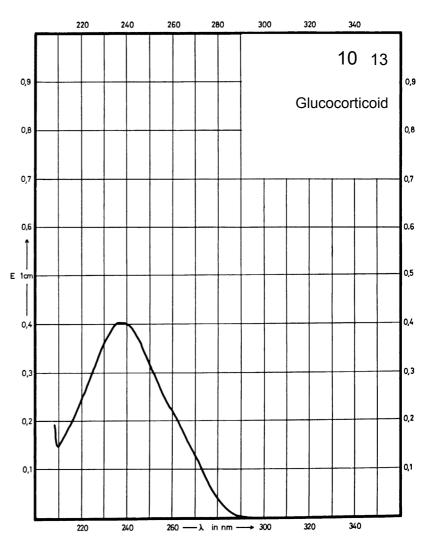

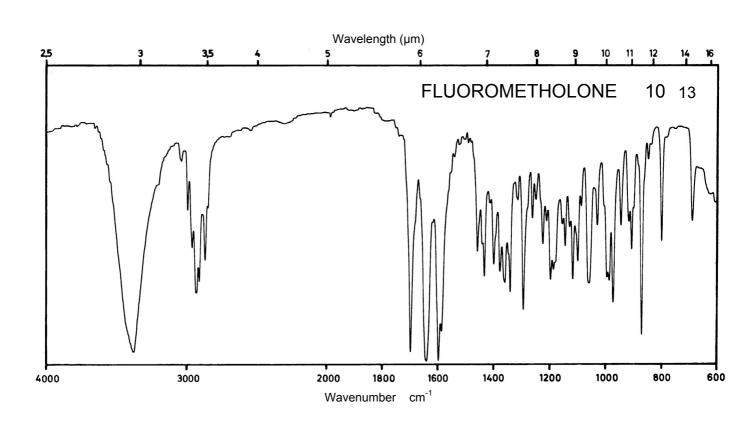




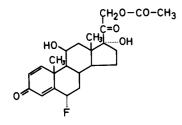


 M_r


420.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm			
E ^{1%} 1cm	381			
ε	16020			

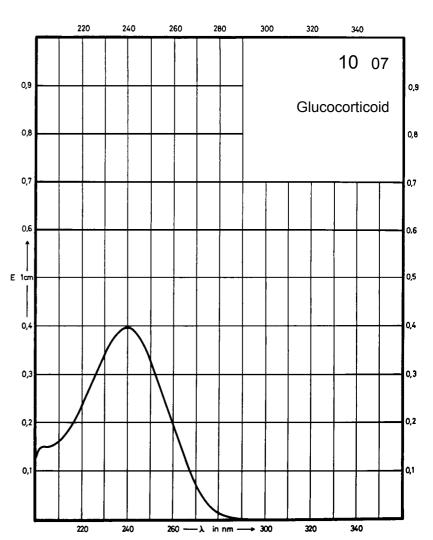

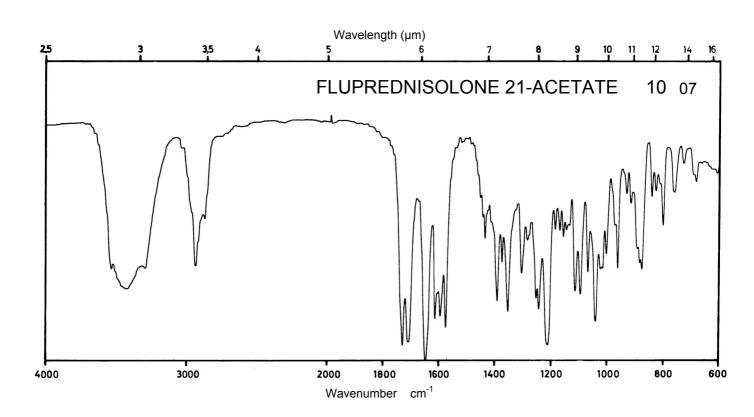





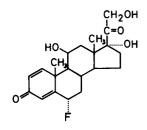

378.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm			
E ^{1%} 1cm	422			
ε	15970			

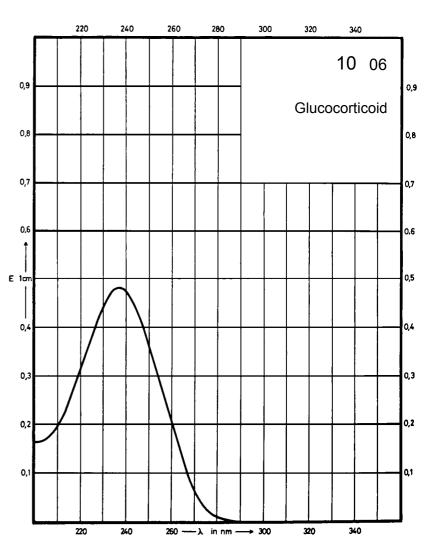

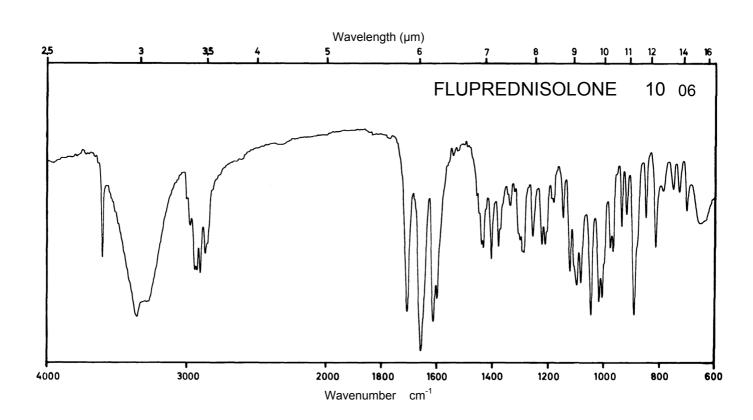







130.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm	265 nm	265 nm	282 nm
E ^{1%} 1cm	537	550	550	420
ε	7000	7200	7200	5500

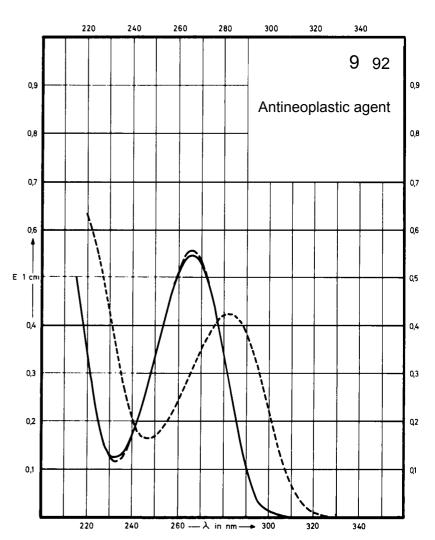

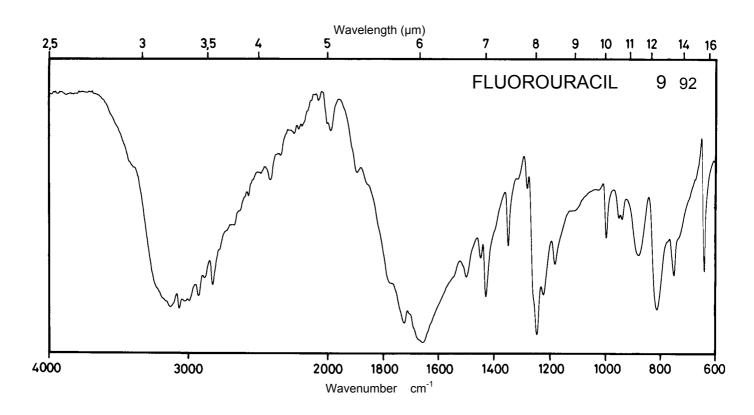







 M_r


129.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 237 nm		285 nm	292 nm
E ^{1%} 1cm	403 560		710	608
ε	5200 7200		9200	7850

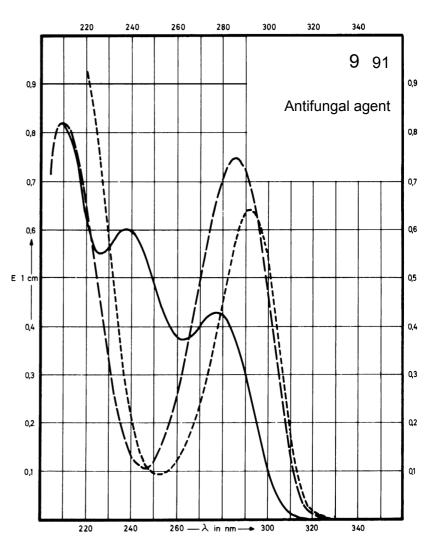

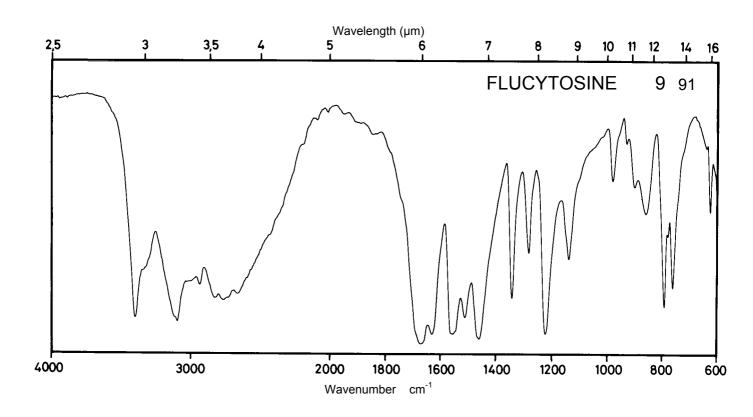





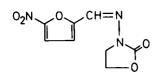


225.2

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of absorption	355 nm 260 nm		366 nm 259 nm	Decom- position observed
E ^{1%} 1cm	840 450		386 737	
ε	18920 10130		16890 12950	

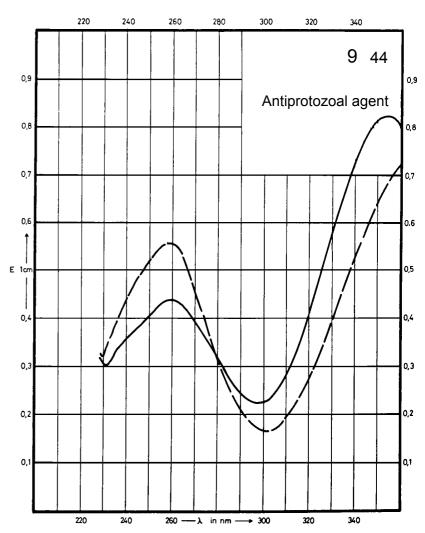

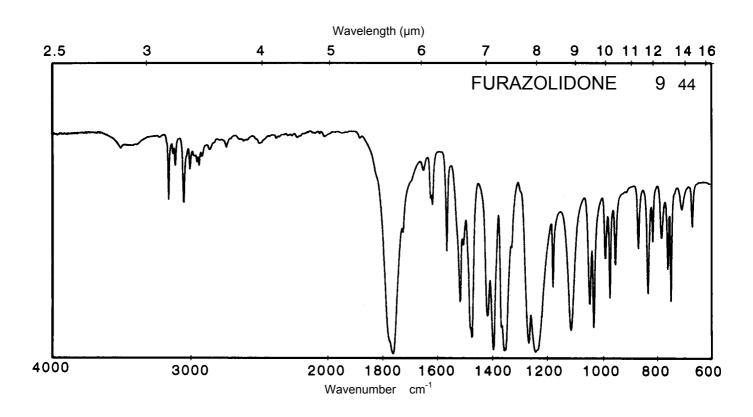





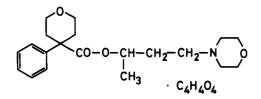


 M_{r} 463.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 259 nm		258 nm	
E ^{1%} 1cm	12.7 17.1		10.7	
3	590 794		496	

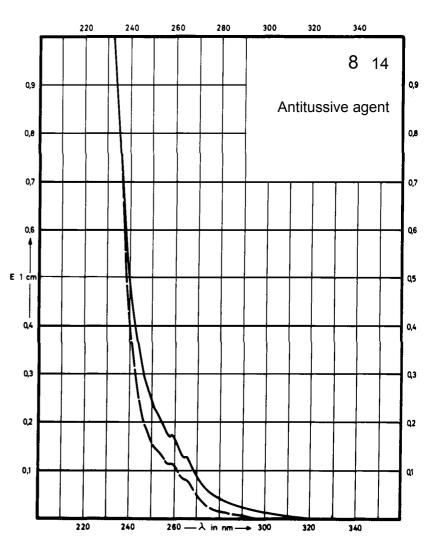

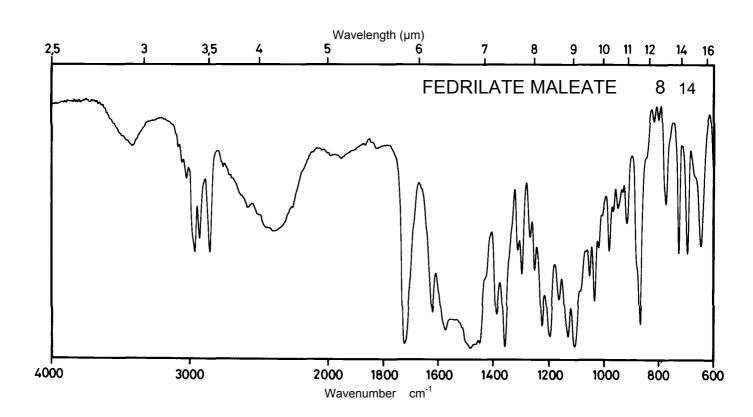



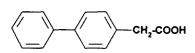



212.3

Concentration


0.5 mg / 100 ml

Solve Symb	 Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximu absorp	252 nm	254 nm	252 nm	254 nm
E ^{1%}	998	1005	969	1002
3	21200	21300	20600	21300

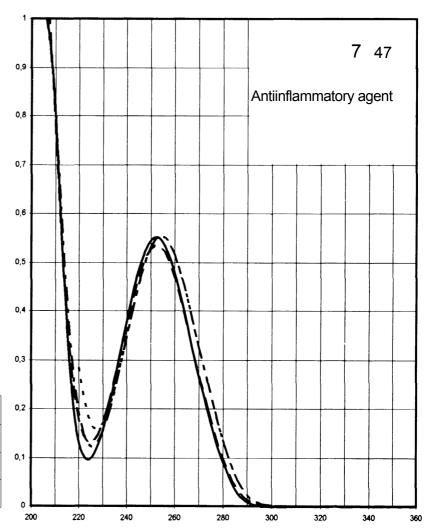

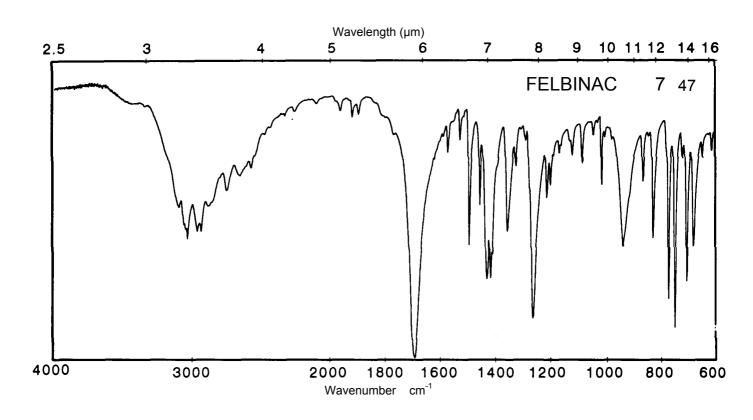





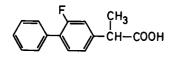


 M_r 244.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	247 nm	247 nm	247 nm	247 nm
E ^{1%} 1cm	808	798	787	804
ε	19700	19500	19200	19600

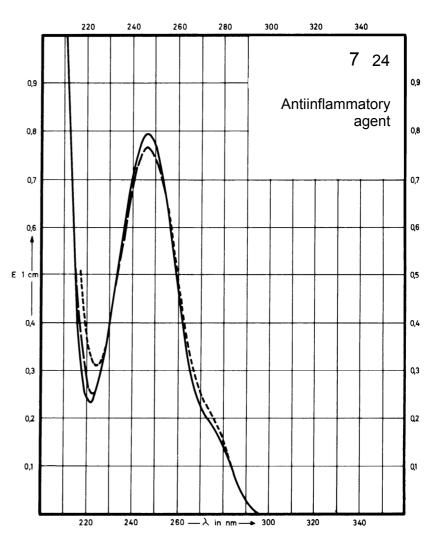

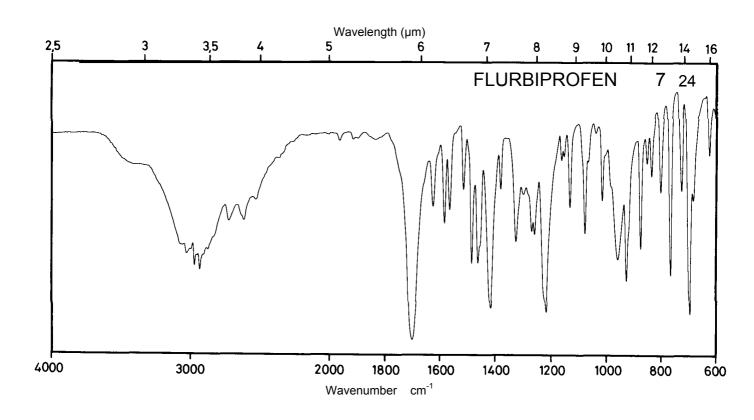





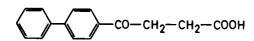


 M_r 254.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm			286 nm
E ^{1%} 1cm	833			816
3	21170			20740

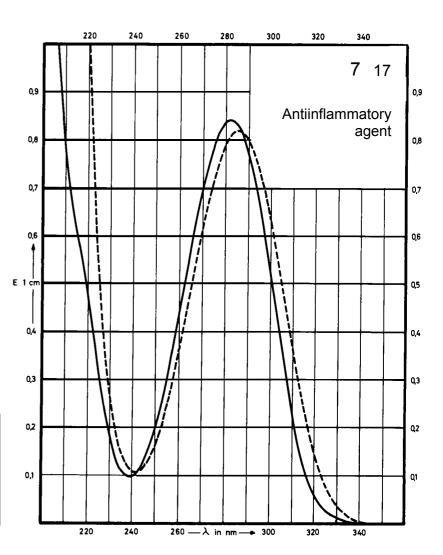

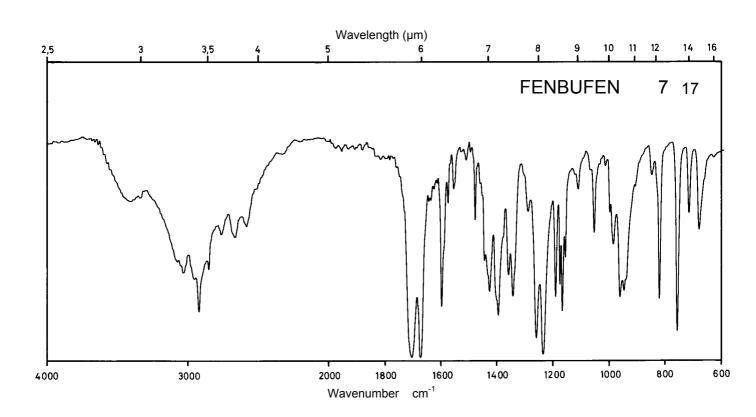


Name FENDOSAL

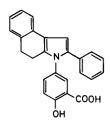


 M_{r}


381.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	292 nm			292 nm
E 1% 1cm	566			564
ε	21590			21510

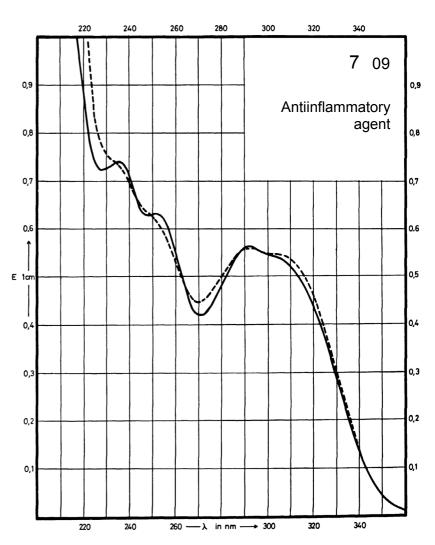

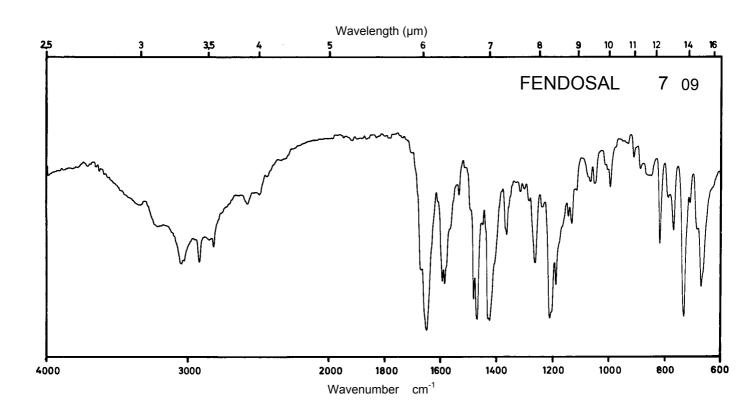





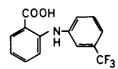


281.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	340 nm 287 nm			287 nm
E ^{1%} 1cm	277 587			530
ε	7790 16510			14900

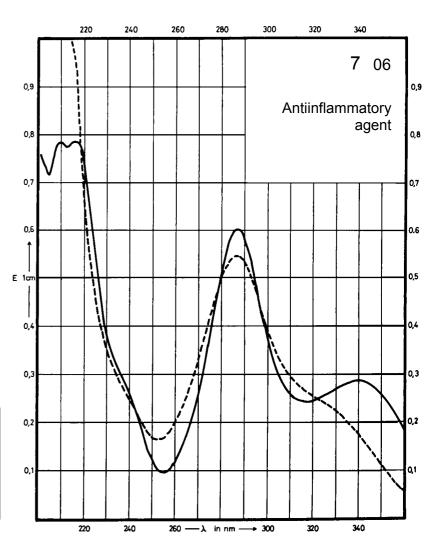

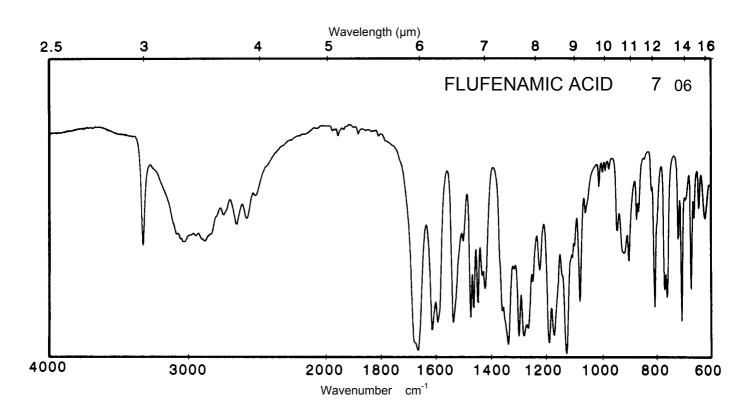


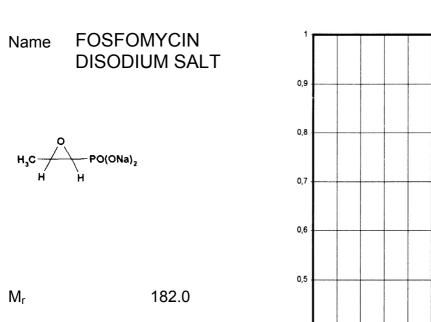


50 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

Concentration

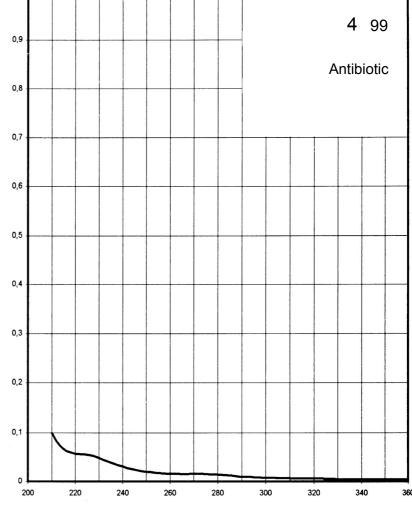

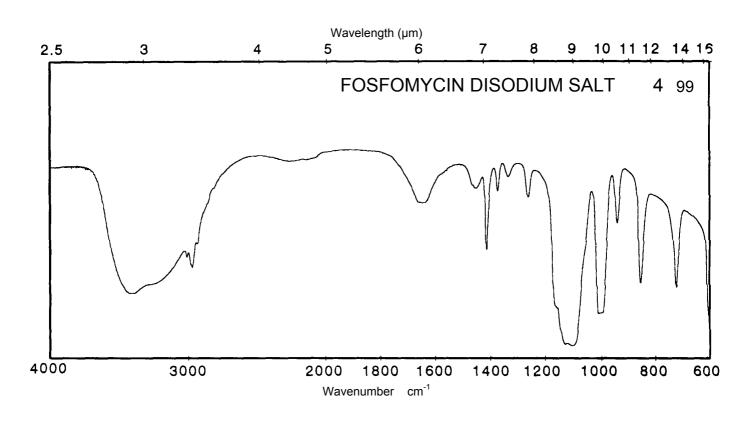





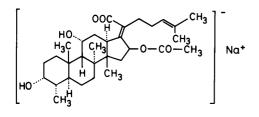


538.7

_					
	Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
	Maximum of absorption				
	E ^{1%} 1cm				
Ī	3				

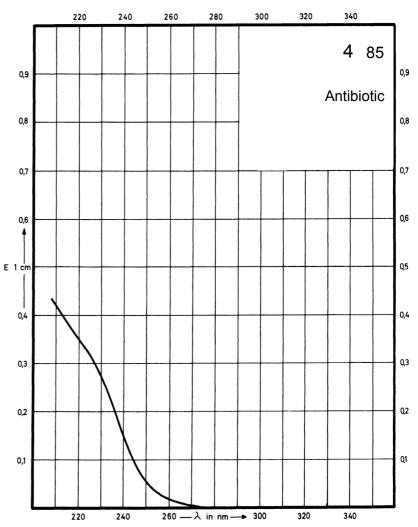

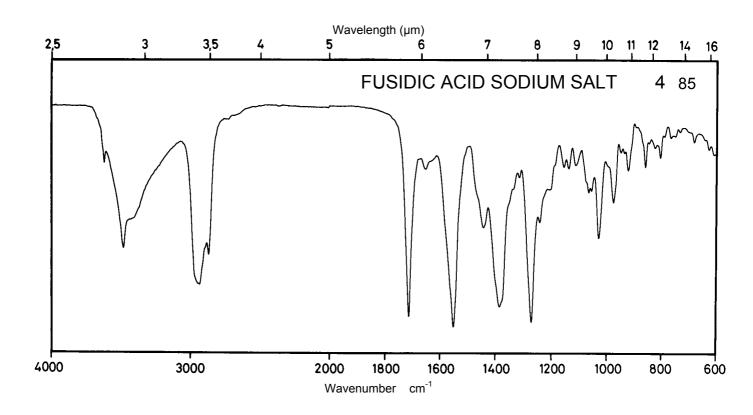




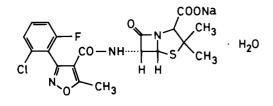


 M_r 493.9

Concentration 1 mg / 100 ml 26 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm 268 nm		Decom- position observed	
E ^{1%} 1cm	21 22			
ε	1050 1075			

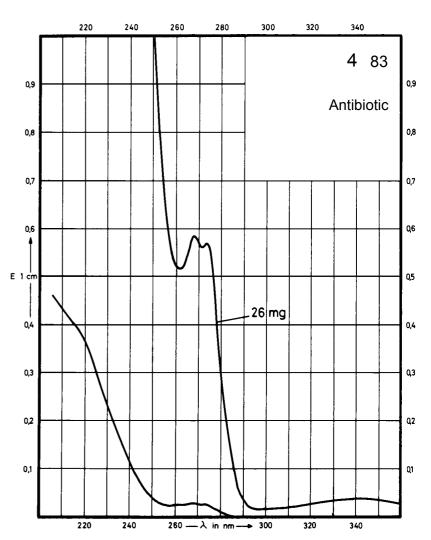

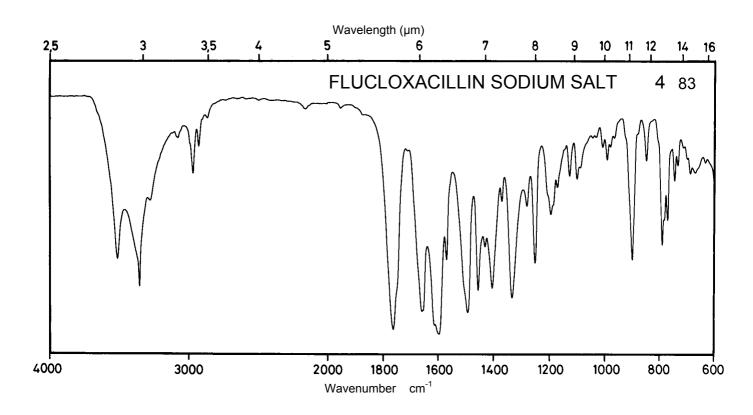





C19H20CIN3 · C16H18N2O4S

 M_r


660.2

Ī	Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
ŀ		283 nm	281 nm		
	Maximum of absorption	275 nm 252 nm	274 nm 252 nm		
ŀ	10/	94	88		
	E ^{1%} 1cm	96	107		
L	ICIII	128	122		
		6200	5800		
	3	6360	7040		
	-	8470	8050		

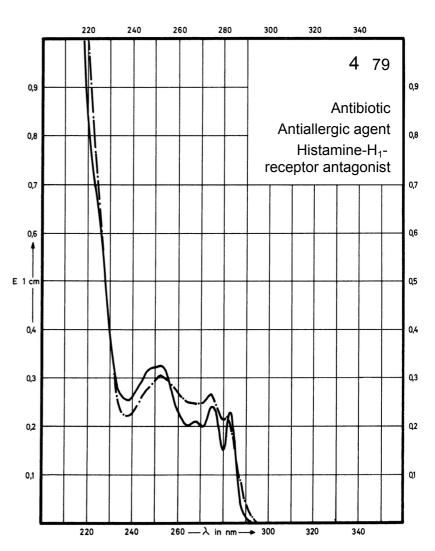

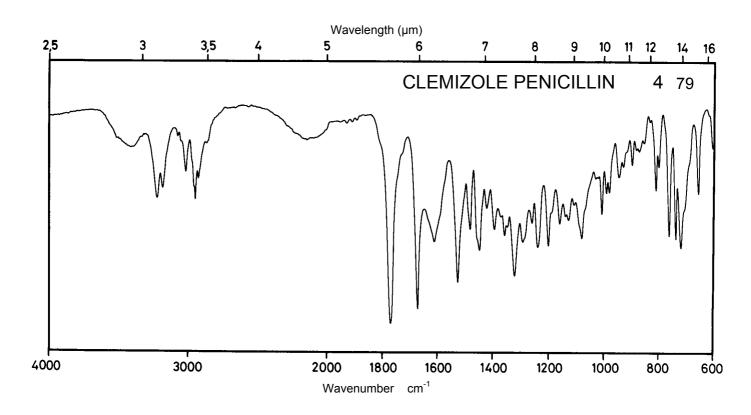


Name FEXOFENADINE HYDROCHLORIDE

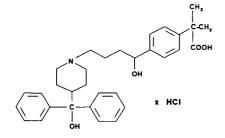


 M_{r}


501.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

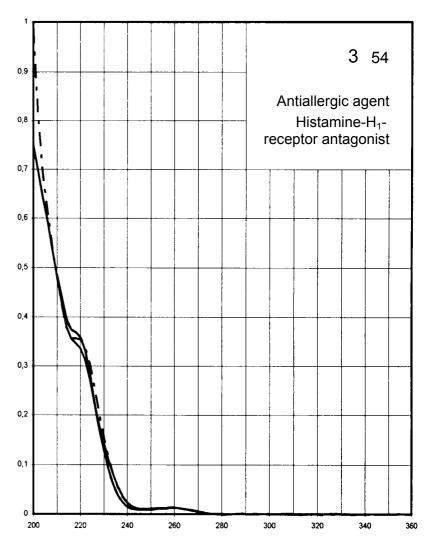

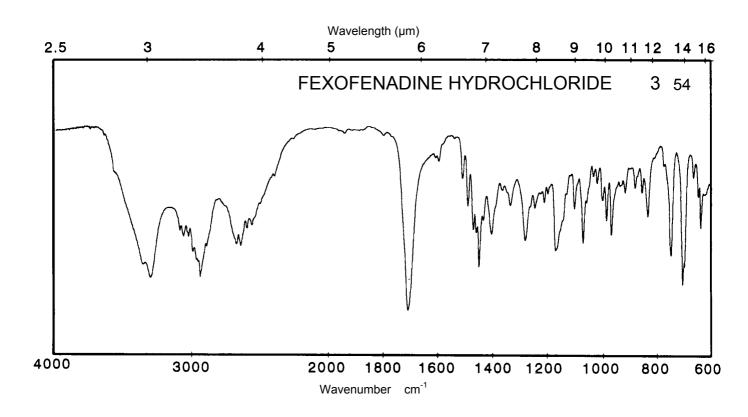



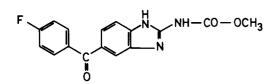



313.3

Concentration


0.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			283 nm 232 nm	
E ^{1%} 1cm			566 850	
ε			17700 26600	

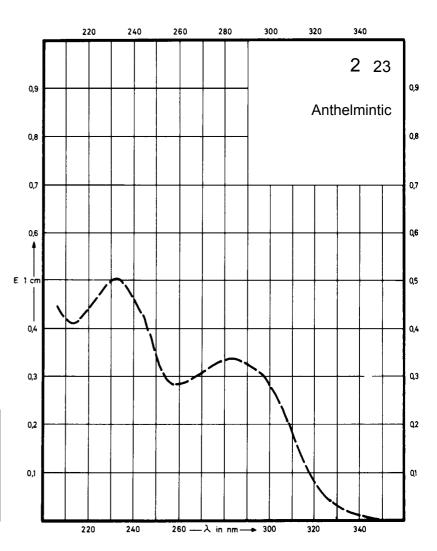

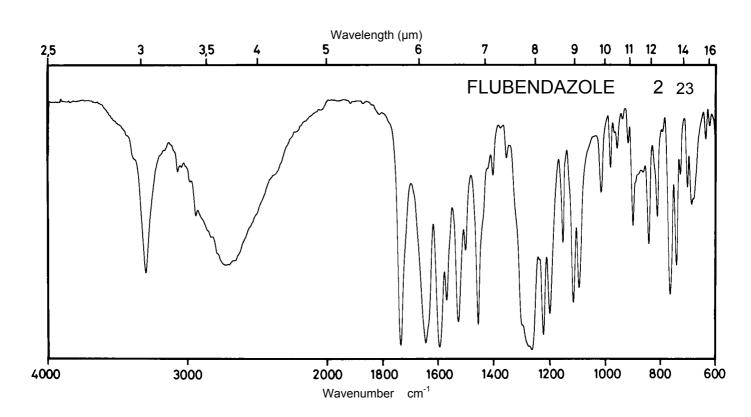





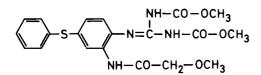

 M_r 446.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm	278 nm	280 nm	282 nm 246 nm
E ^{1%} 1cm	353	305	256	388 650
ε	15700	13600	11400	17300 29000

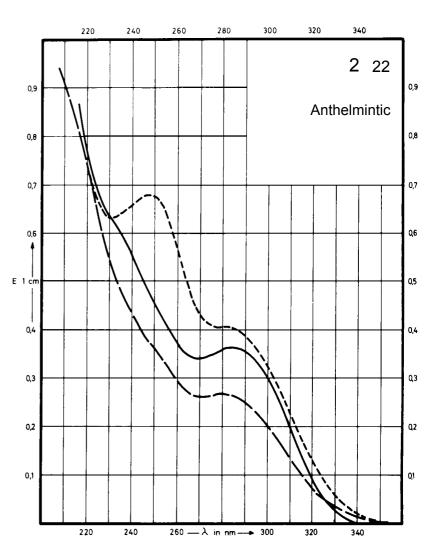


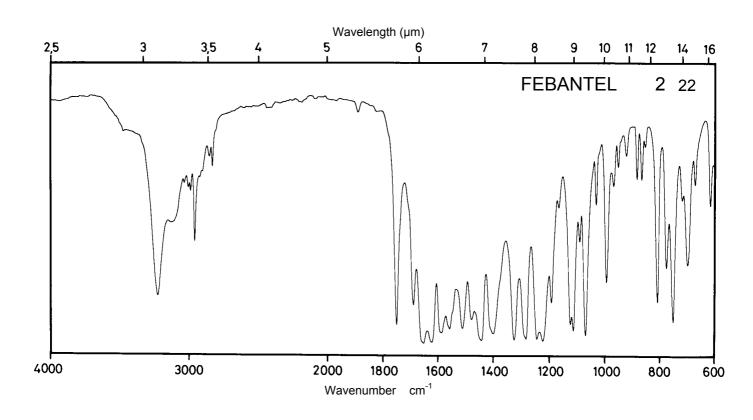


Name FENBENDAZOLE


 M_{r} 299.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	294 nm		290 nm	
E ^{1%} 1cm	490		390	
ε	14680		11680	

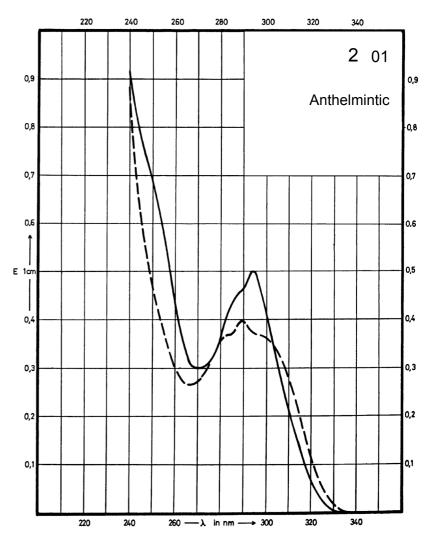

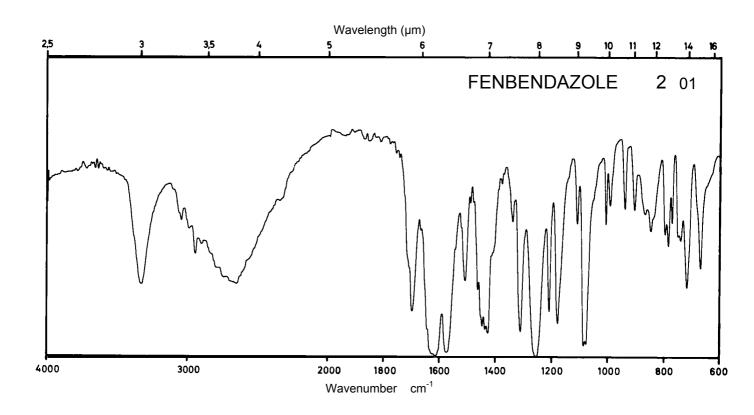




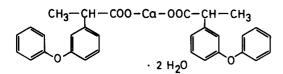


 M_r 558.6


Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm 273 nm		272 nm	272 nm
E ^{1%} 1cm	61 69		62	63
ε	3420 3840		3460	3520

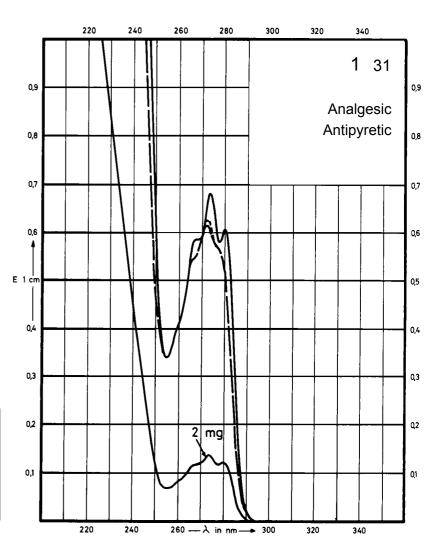

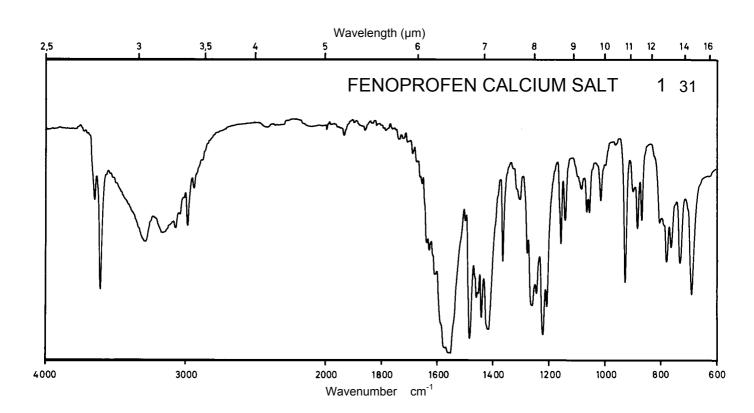


Name FENTANYL CITRATE

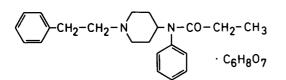


 M_{r}


528.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 258 nm 252 nm	263 nm 258 nm 252 nm	263 nm 258 nm 252 nm	
E ^{1%} 1cm	5.98 8.06 7.58	6.42 8.33 7.34	6.42 8.33 7.34	
ε	316 426 400	340 440 390	340 440 390	

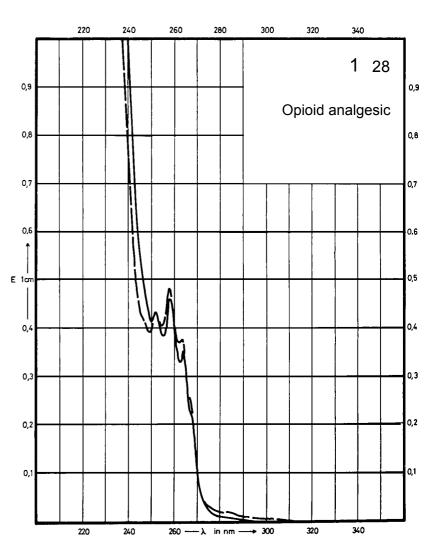

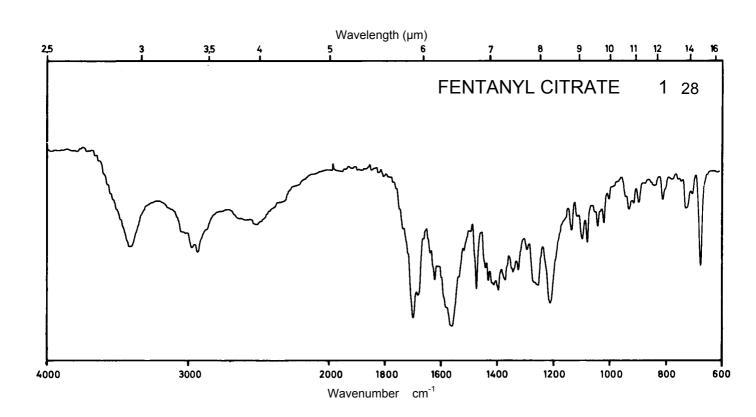


Name FLOCTAFENINE

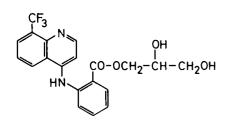


 M_{r}


406.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	362 nm		348 nm 233 nm	359 nm 249 nm
E 1% 1cm	477		451 545	390 392
ε	19380		18320 22150	15860 15950

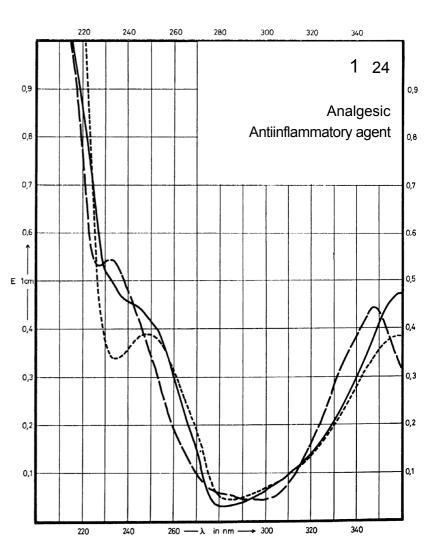

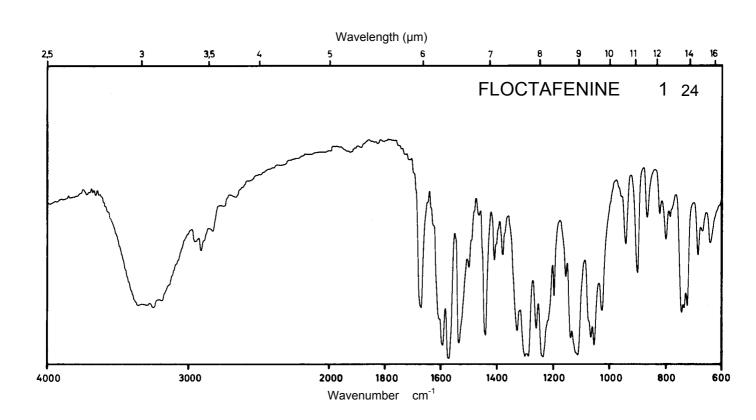





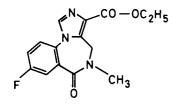


 M_r 303.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm			Decom- position observed
E 1% 1cm	684			
ε	20740			

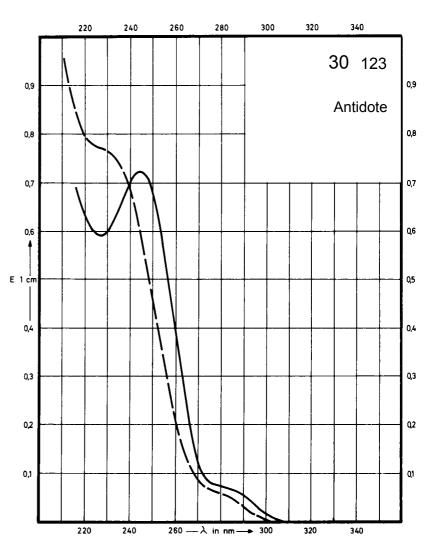

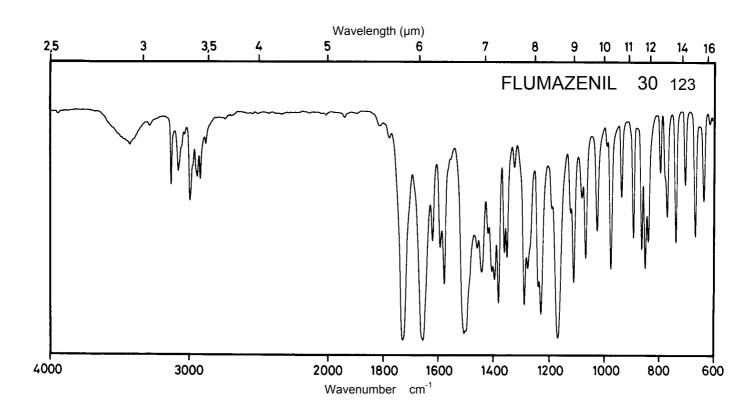





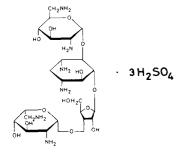

908.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

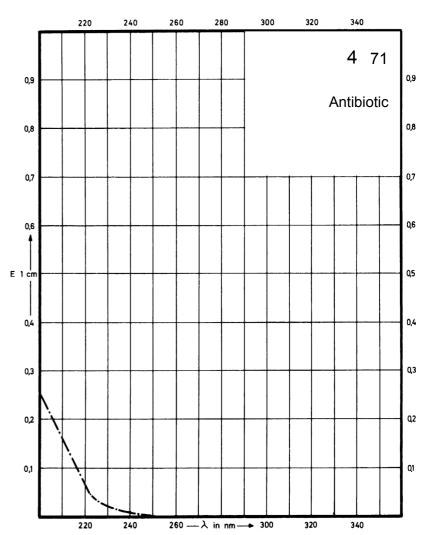

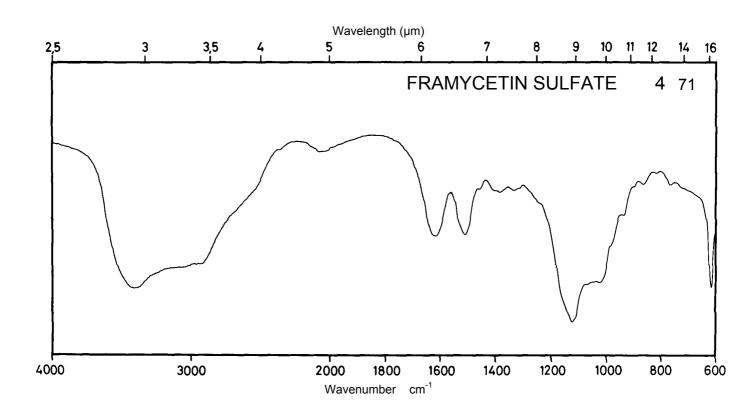







154.1

Concentration 1 m


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	333 nm		328 nm	
E ^{1%} 1cm	271		261	
ε	4180		4020	

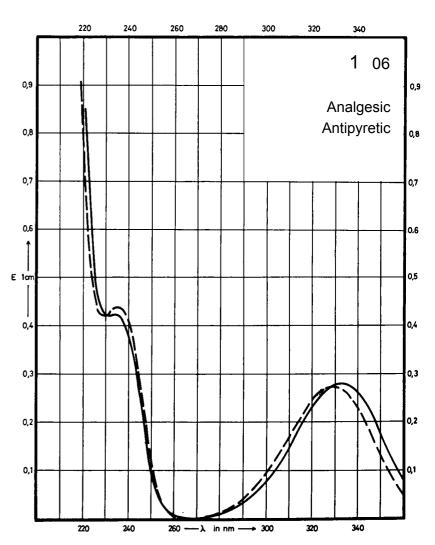

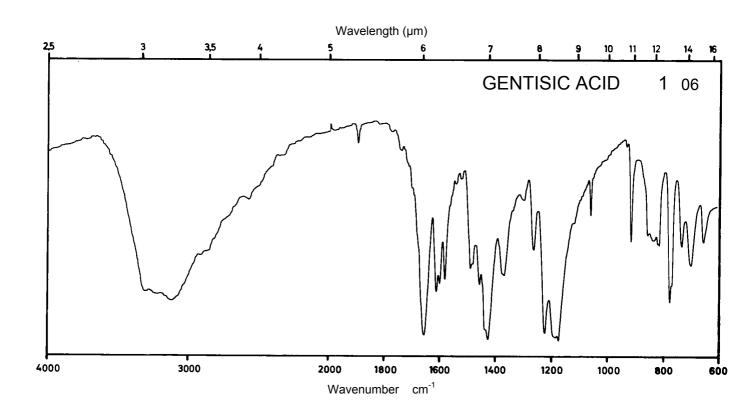





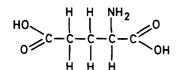


147.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

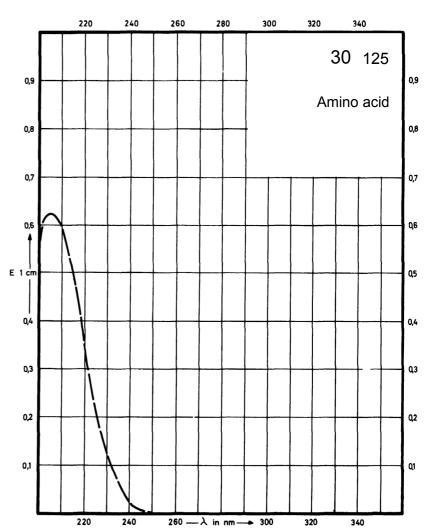

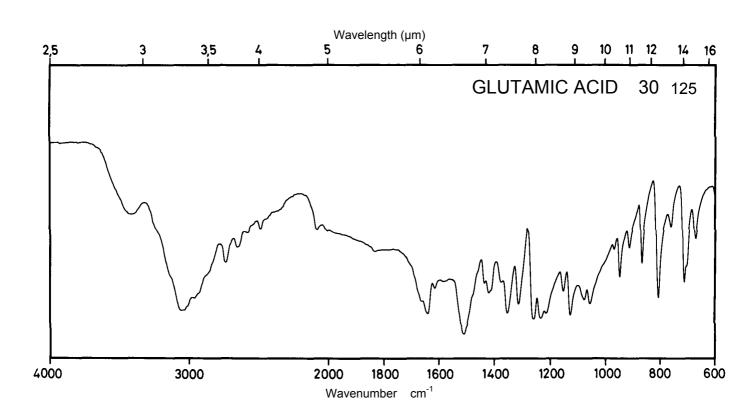



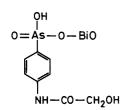



499.1

Concentration


1.25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			254 nm	
E ^{1%} 1cm			423	
3			21100	

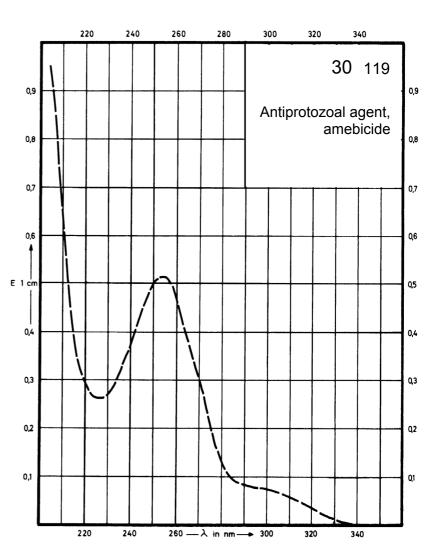

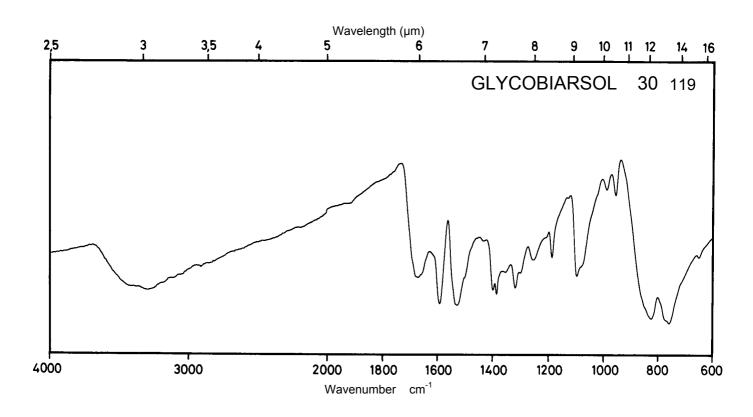



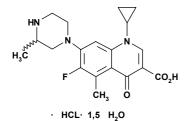


 M_r


422.9

Concentration 0


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	319 nm 280 nm	321 nm 278 nm 218 nm	320 nm 279 nm 218 nm	326 nm 273 nm
E ^{1%} 1cm	326 1060	294 975 367	295 1010 385	313 869
ε	13800 44800	12400 41200 15500	12500 42700 16300	13300 36700

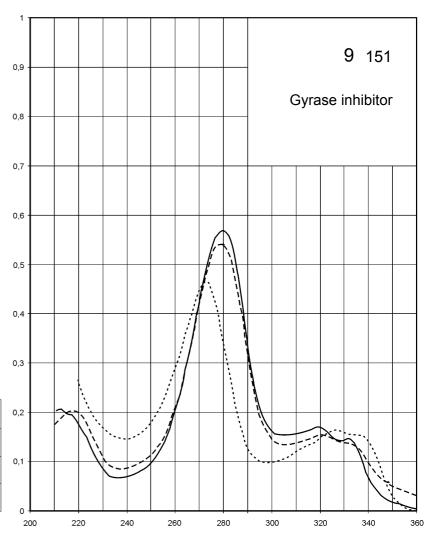

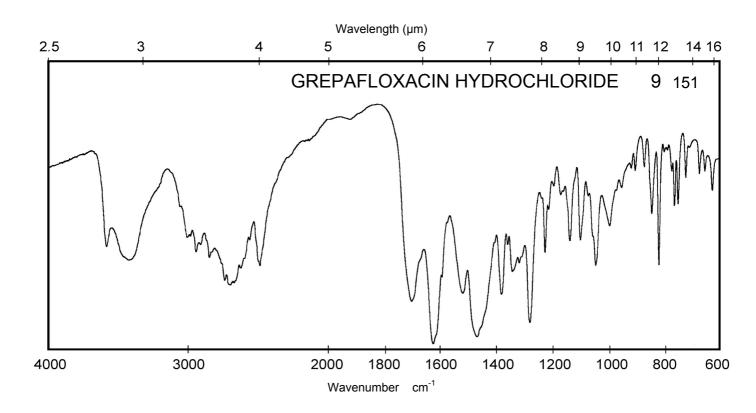




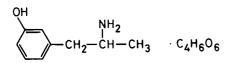


 M_r 301.3


Concentration 3 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm		272 nm	289 nm 237 nm
E ^{1%} 1cm	69		57	97 291
ε	2100		1720	2900 8800

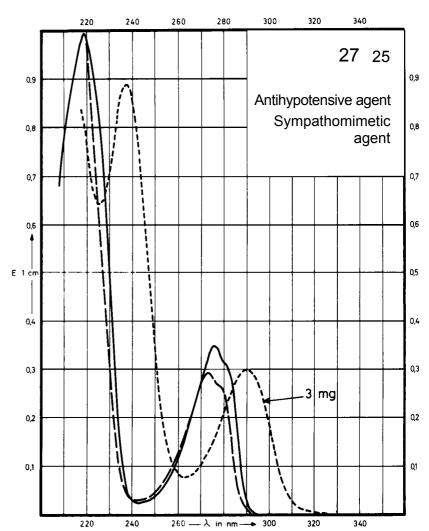

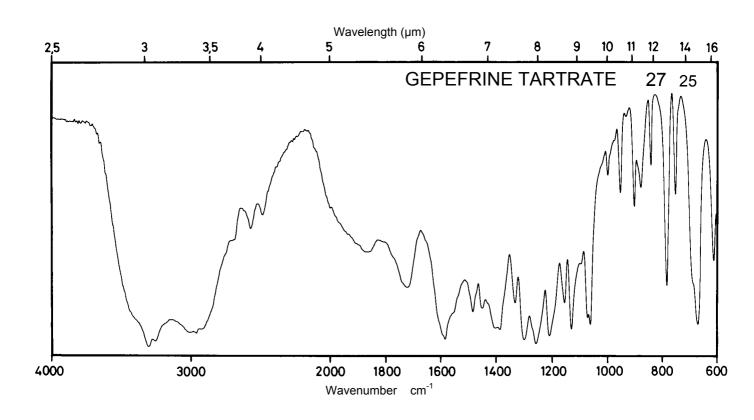




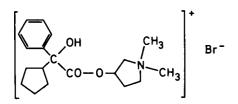


398.4

Concentration 10 mg / 100 ml 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm	263 nm 257 nm 251 nm		
E ^{1%} 1cm	4.36 5.64 5.52	4.83 5.90 5.47		
ε	174 225 220	193 235 218		

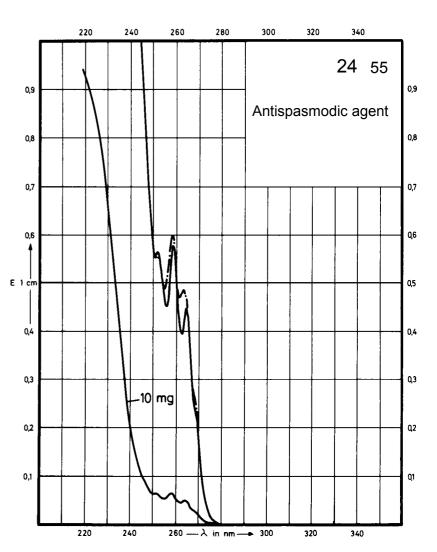

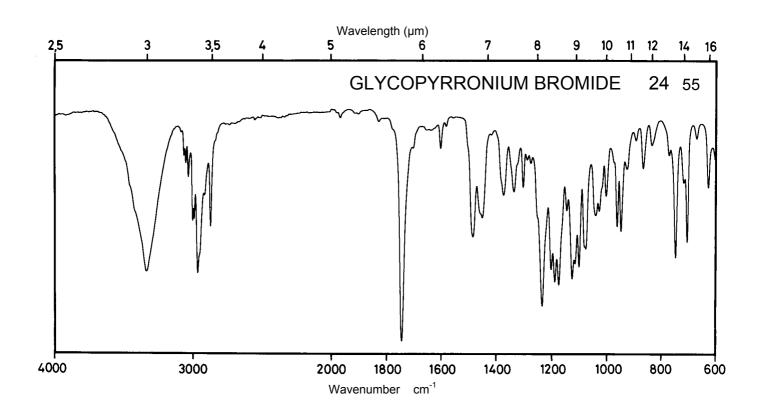

Name GUAIFENESIN

0-сH₂-сH(OH)-сH₂OH 0-сH₃


 M_{r}


198.2

Concentration

2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm 222 nm		273 nm 221 nm	272 nm
E ^{1%} 1cm	122 371		120 370	121
3	2420 7350		2380 7330	2400

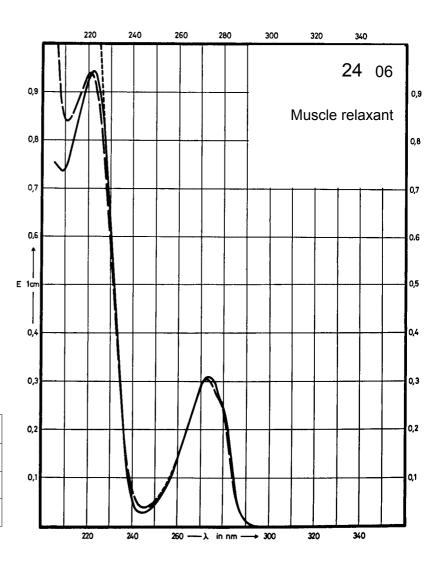

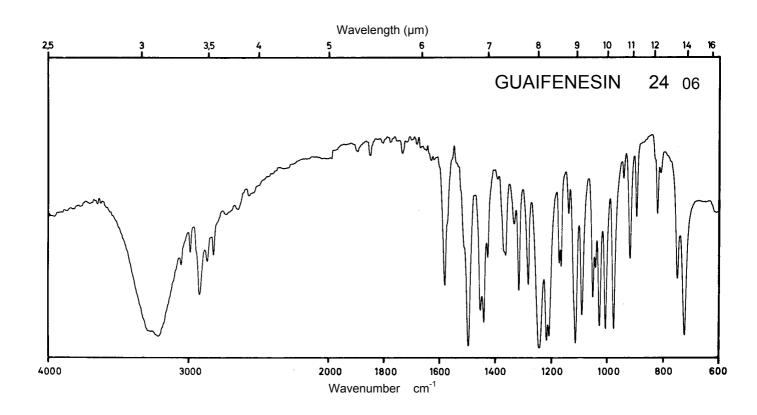

Name GALLAMINE TRIETHIODIDE

$$\begin{bmatrix} O-CH_2-CH_2-N(C_2H_5)_3 \\ O-CH_2-CH_2-N(C_2H_5)_3 \\ O-CH_2-CH_2-N(C_2H_5)_3 \end{bmatrix}^{+++}$$


 M_r


891.5

Concentration

60 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			266 nm 223 nm	265 nm
E ^{1%} _{1cm}			8.4 518	8.4
ε			750 46180	750

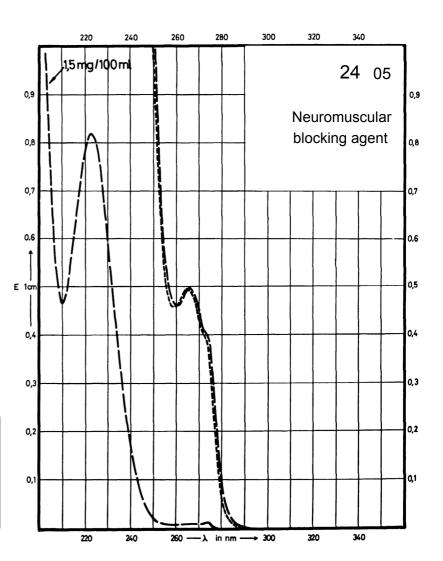

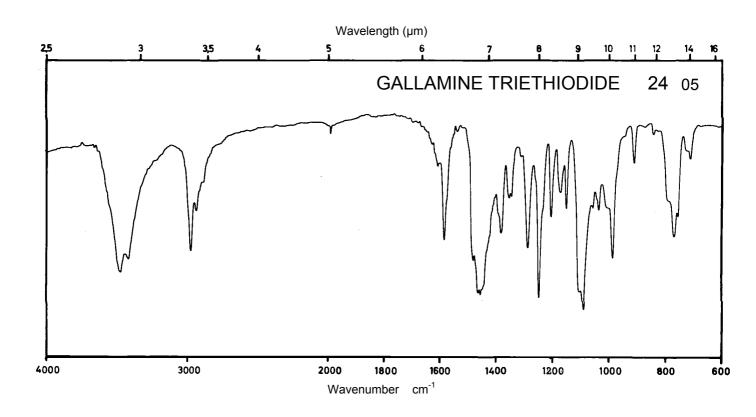



O-C(CH₃)₃


|
Glp-His-Trp-Ser-Tyr- p-Ser-Leu-Arg-Pro
|
H₂N-CO-NH-NH


x CH3COOH

 M_r 1329.5

Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm 280 nm 220 nm	279 nm 218 nm	278 nm 218 nm	289 nm 281 nm
E ^{1%} 1cm	40.6 53.1 374	45.0 336	47.4 340	47.7 49.6
ε	5400 7050 49700	5600 44700	6300 45100	6340 6600

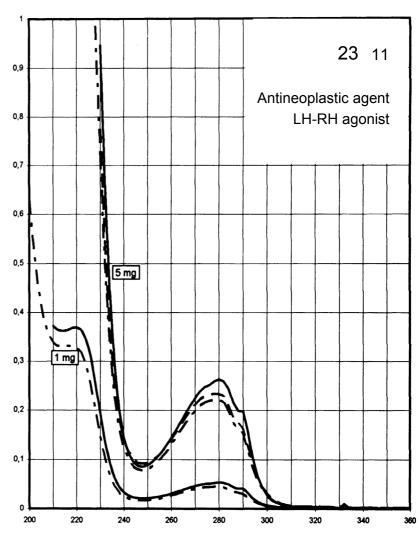

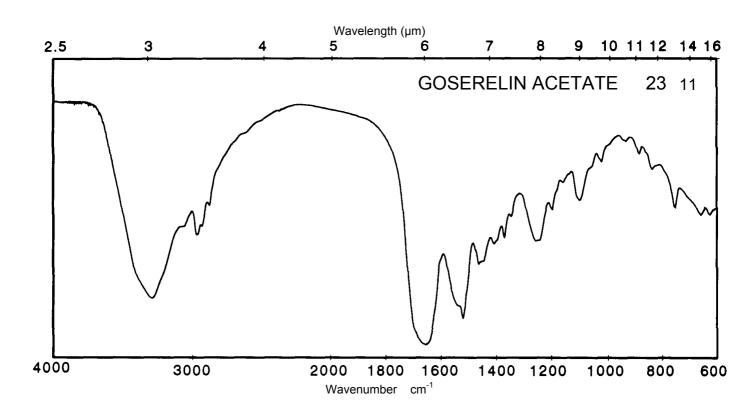



23 06 0,9 Hormone 0,8 0,7 0,7 0,6 0,6 0,5 E 1 cm 0,4 0,4 0,3

280

260

220

0,3

0,2

0,1

300


320

340


 M_{r} 1182.3

10 mg / 100 ml Concentration

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm	279 nm	279 nm	282 nm
E ^{1%} 1cm	54	49	49	52
3	6360	5780	5780	6100

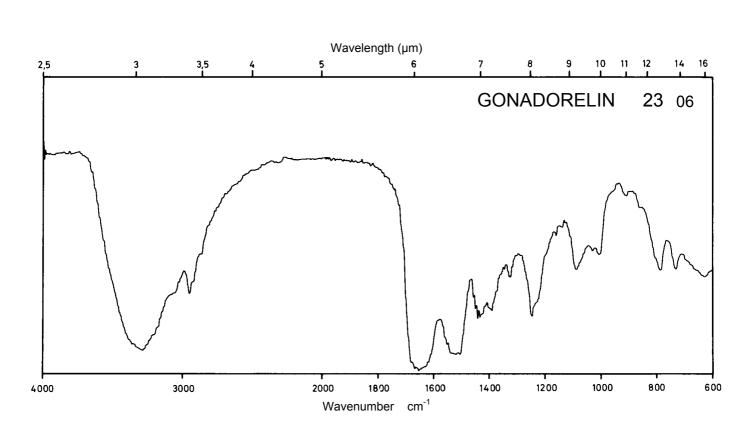





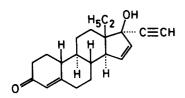


 M_r 310.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm	247 nm		
E ^{1%} 1cm	569	569		
3	17660	17660		

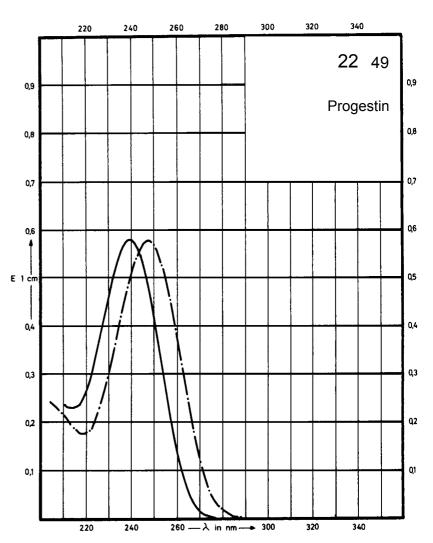

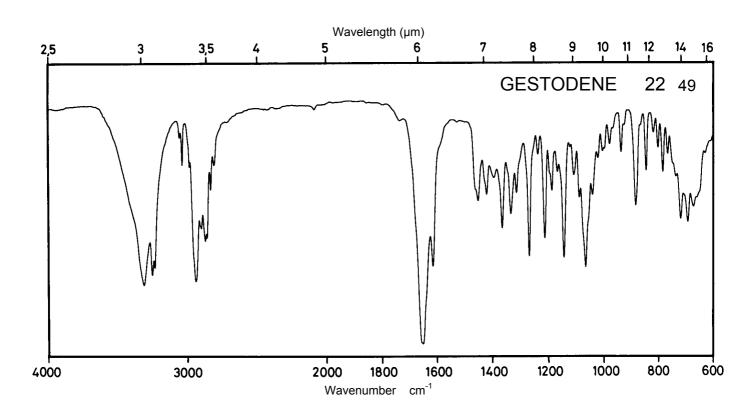



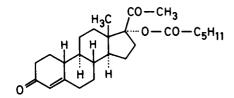



414.6

Concentration


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm			
E ^{1%} 1cm	423			
ε	17540			

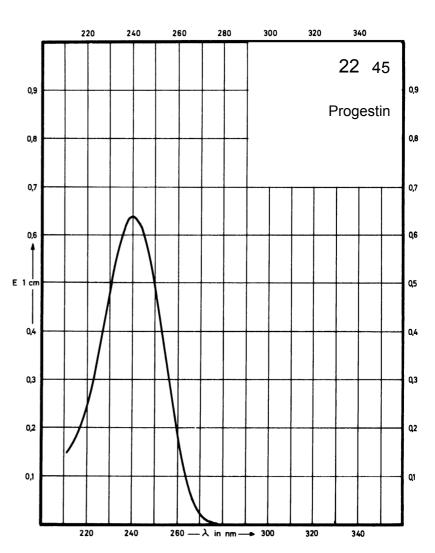

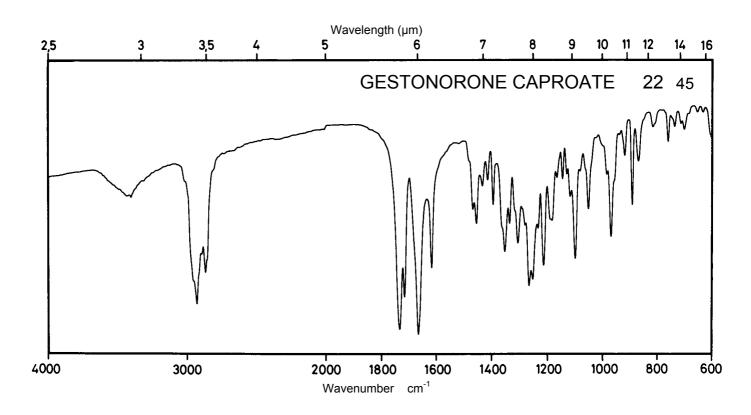





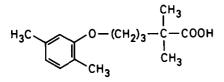

 M_r 250.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 219 nm		273 nm	
E ^{1%} 1cm	80 443		66	
ε	2000 11100		1660	

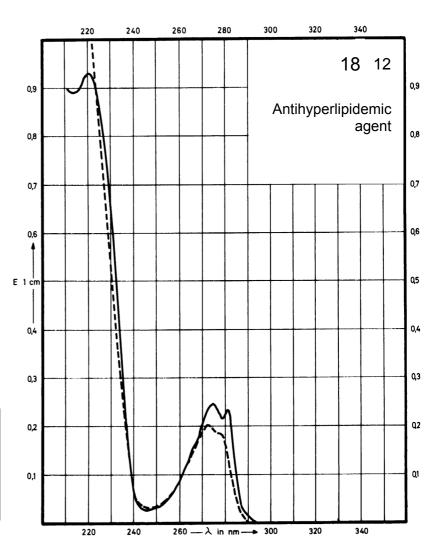

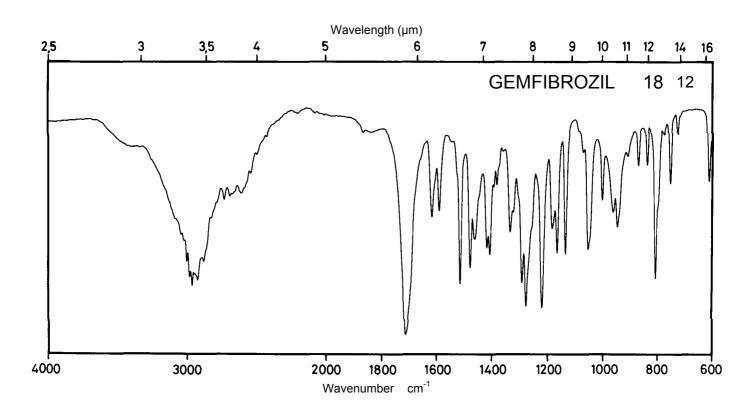


Name GLUTETHIMIDE


0 N O C₂H₅

 M_{r}


217.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	256 nm 250 nm		256 nm 250 nm	258 nm 251 nm
E ^{1%} 1cm	18.2 17.2		19.4 19.7	10.1 9.6
3	395 370		420 430	220 210

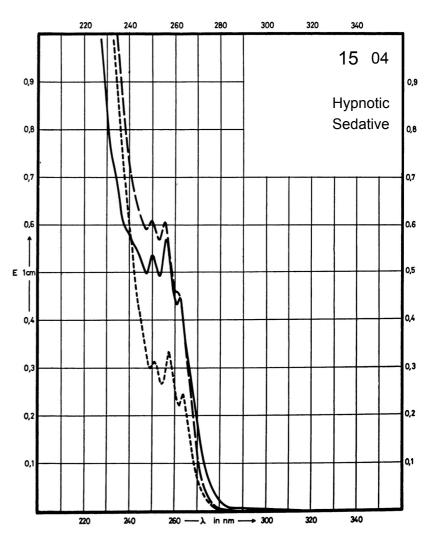

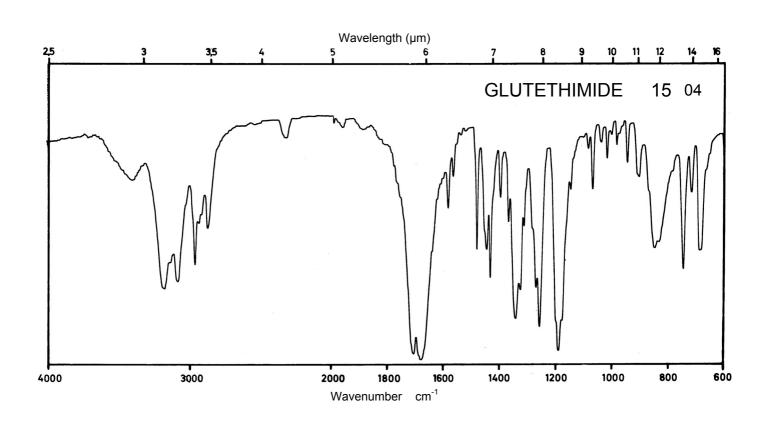




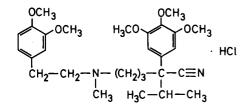


 M_r 521.1


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm		277 nm	277 nm
E ^{1%} 1cm	67		65	63
3	3500		3400	3300

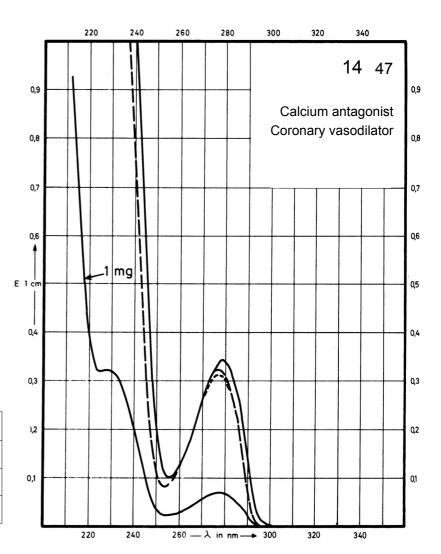

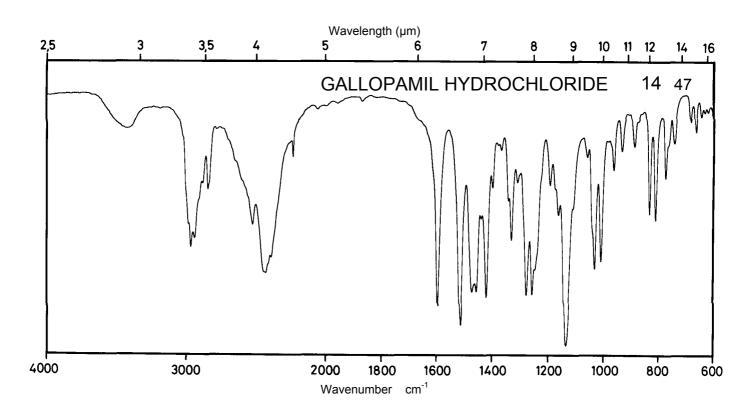

Name GLYCOL SALICYLATE

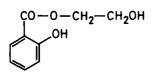


 M_r


182.2

Concentration 1.7


1.7 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	306 nm 237 nm		303 nm 238 nm	332 nm 243 nm
E ^{1%} 1cm	223 480		198 494	Decom- position observed
3	4100 8700		3600 9000	

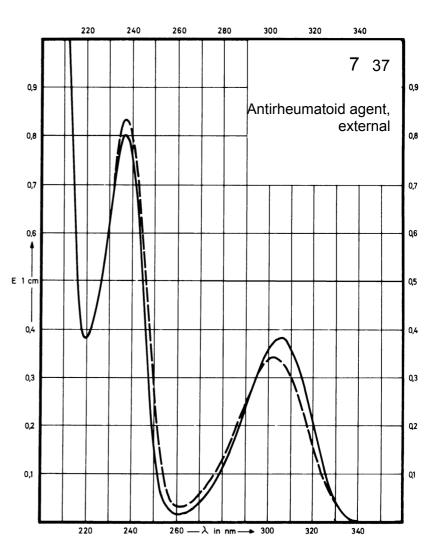

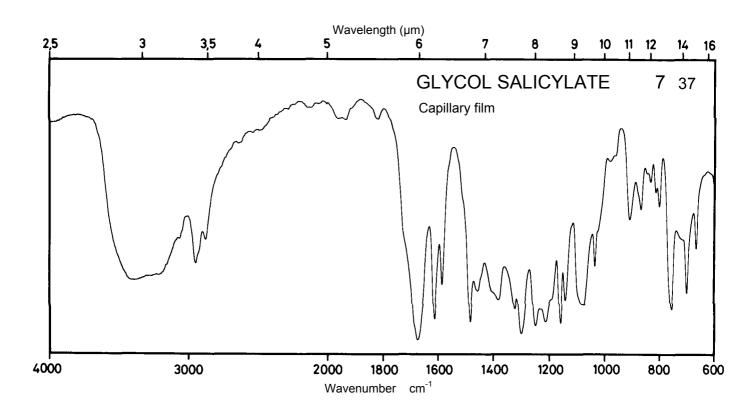





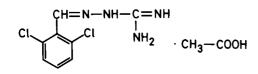


291.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm		271 nm	296 nm
E ^{1%} 1cm	390		455	444
ε	11400		13300	12900

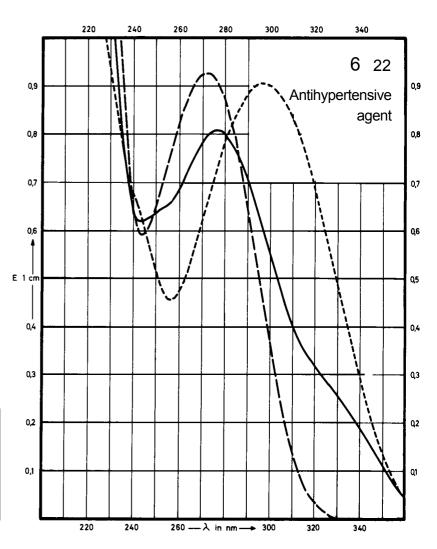

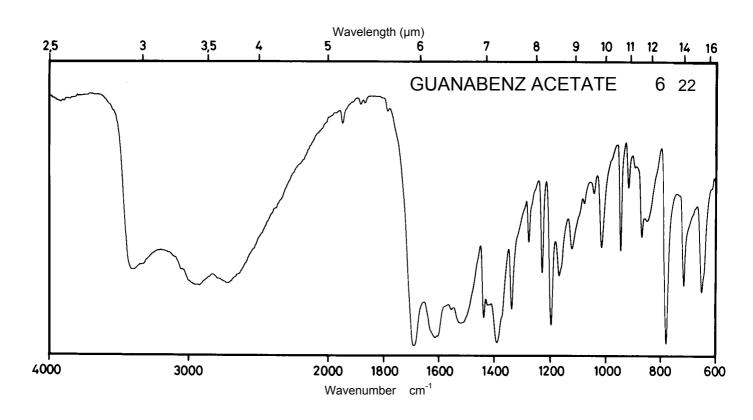





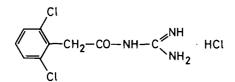


282.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm 272 nm	280 nm 272 nm	280 nm 272 nm	
E ^{1%} 1cm	6.8 8.8	7.8 9.9	7.8 9.8	
ε	192 249	220 280	220 277	

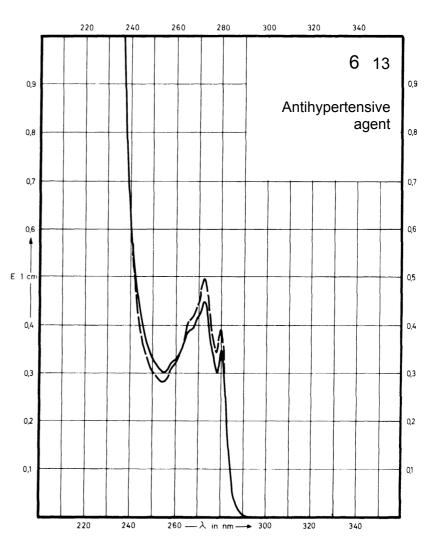

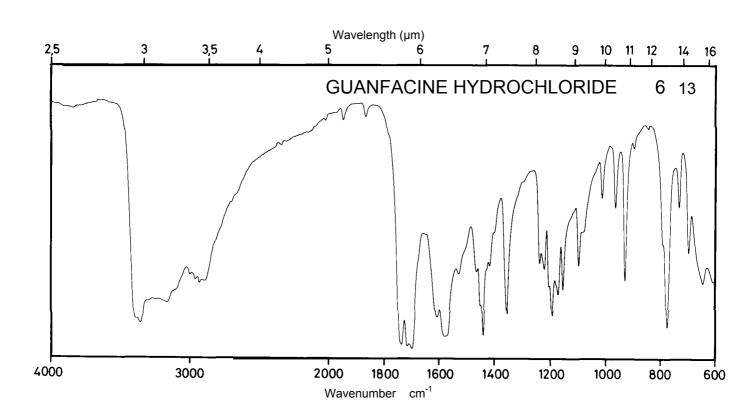





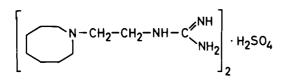


 M_r 494.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

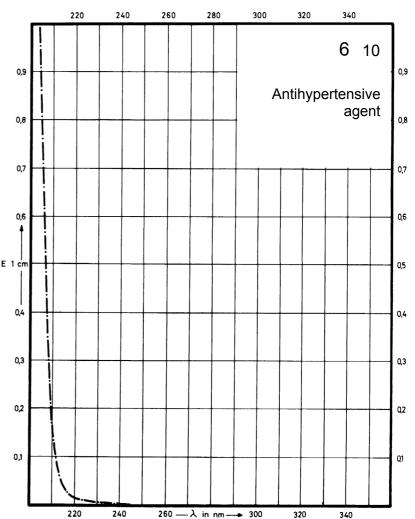

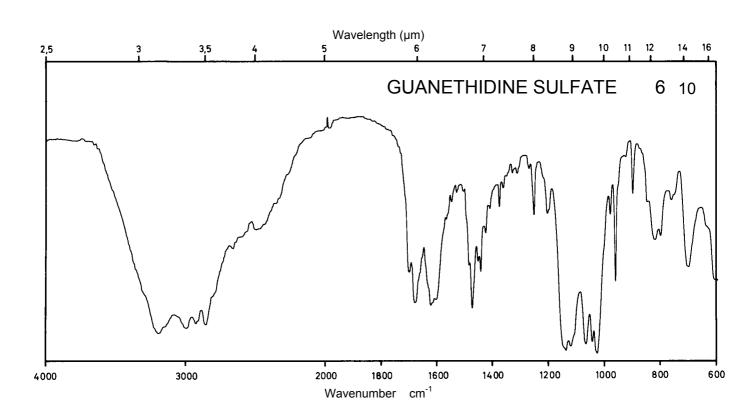





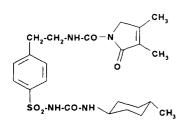


490.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	227 nm			
E ^{1%} 1cm	629			
ε	30800			

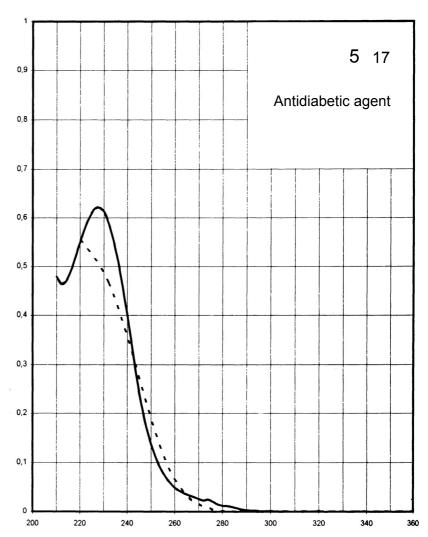

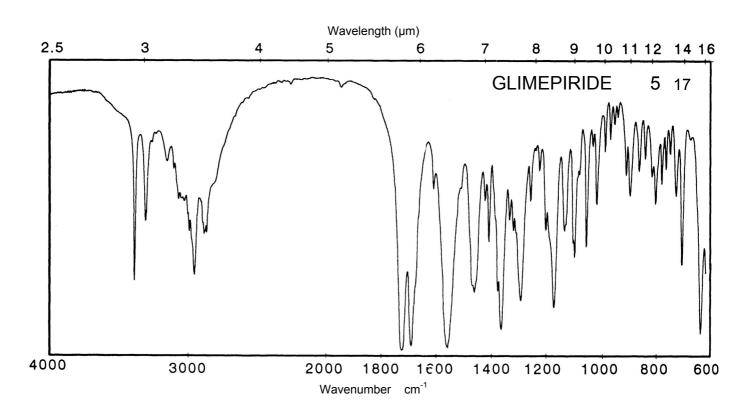

Name GLISOXEPIDE

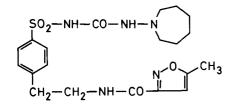



 M_r

449.5


Concentration 1.5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	(274 nm) 229 nm		229 nm	
E ^{1%} 1cm	498		468	
3	22400		21000	

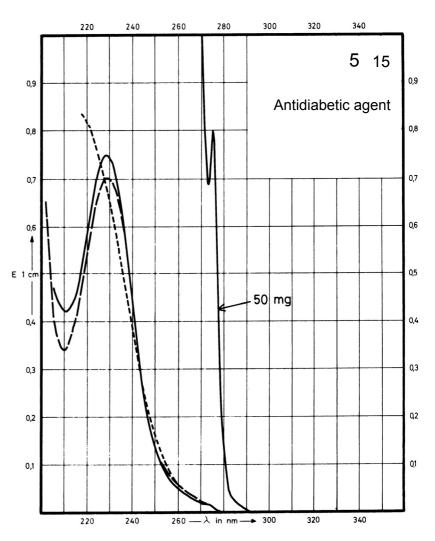

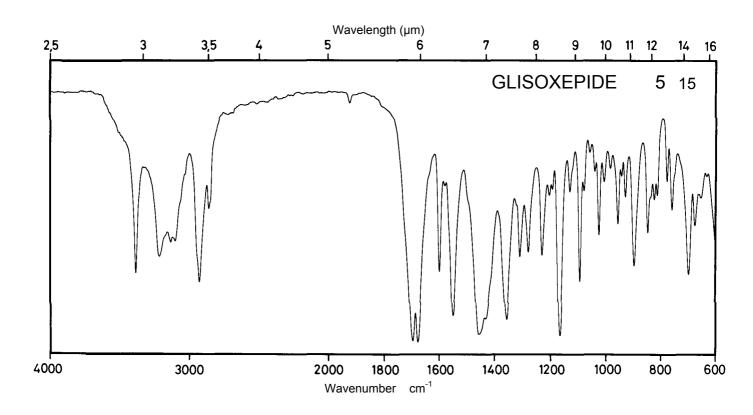





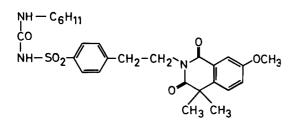


 M_r 527.6

Concentration 1 mg / 100 ml 10 mg / 100 ml

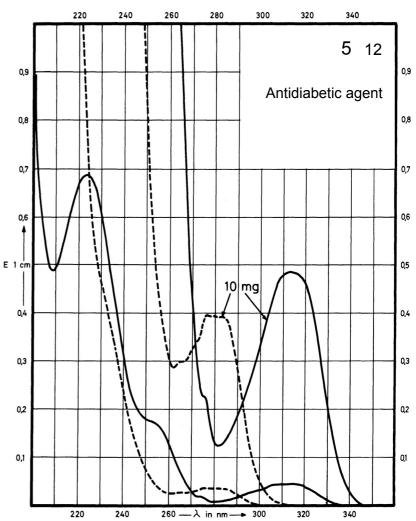

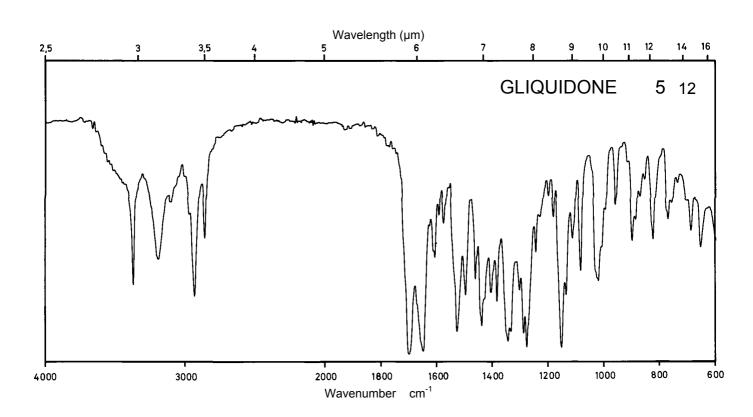




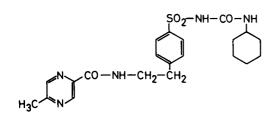


 M_r


445.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 227 nm		276 nm 227 nm	277 nm
E ^{1%} 1cm	239 546		231 513	229
ε	10630 24310		10280 22870	10190

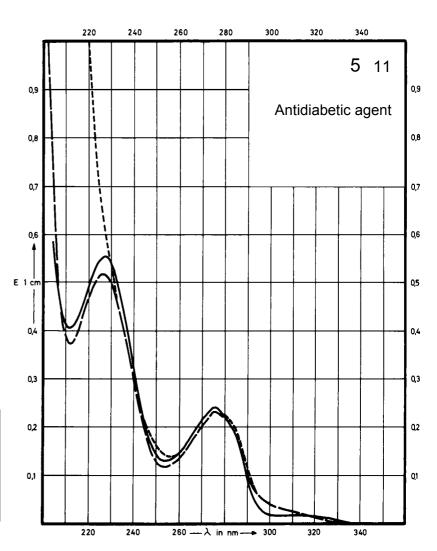

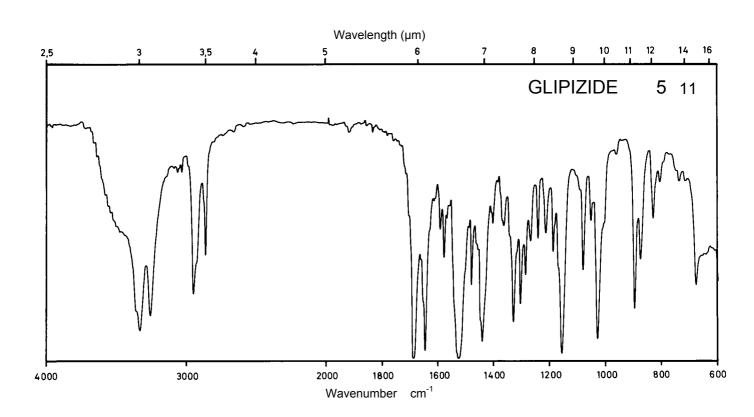




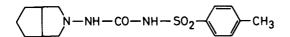


323.4


1 mg / 100 ml Concentration 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 228 nm		230 nm	263 nm
E ^{1%} 1cm	19.7 434		440	19.4
ε	635 14040		14230	625

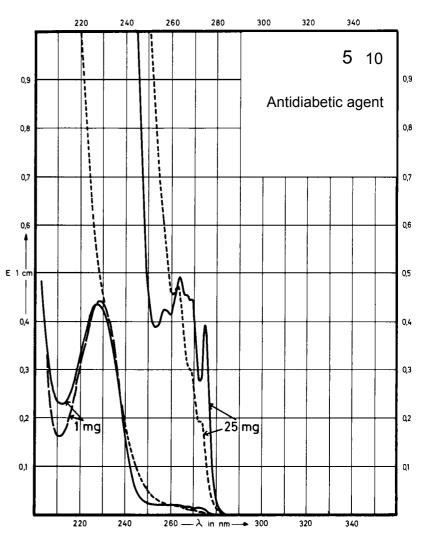

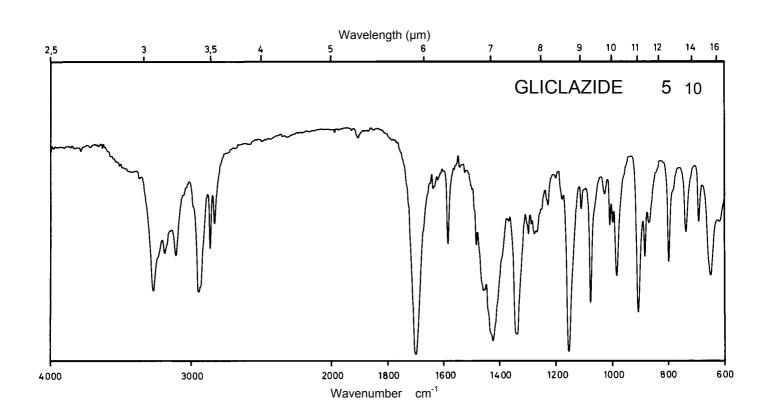





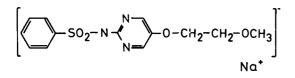


 M_r 331.3

Concentration 1 mg / 100 ml 2 mg / 100 ml

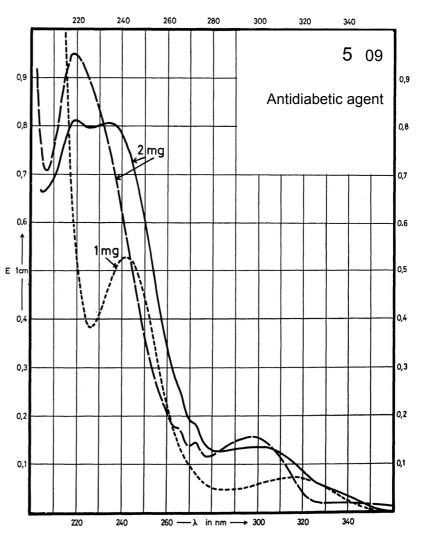

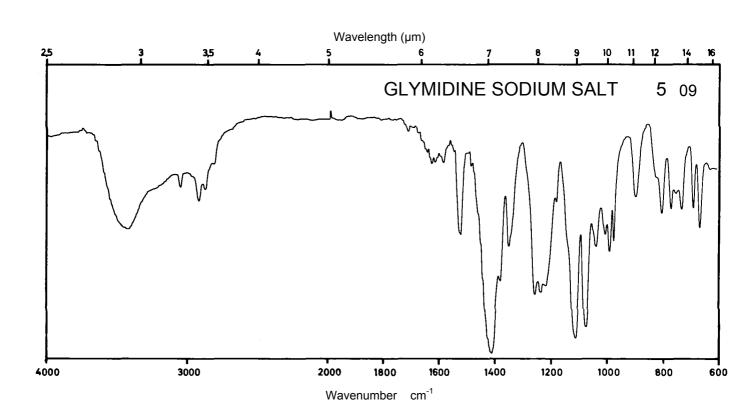





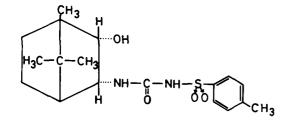



366.5


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm			
E ^{1%} 1cm	16.6			
3	610			

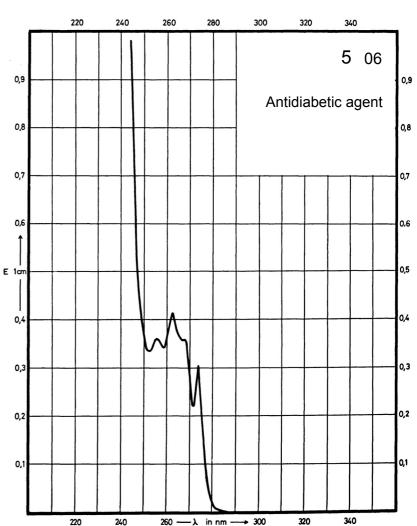


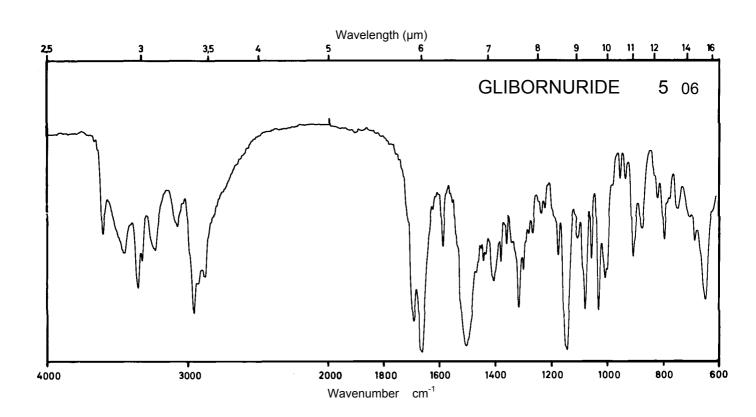


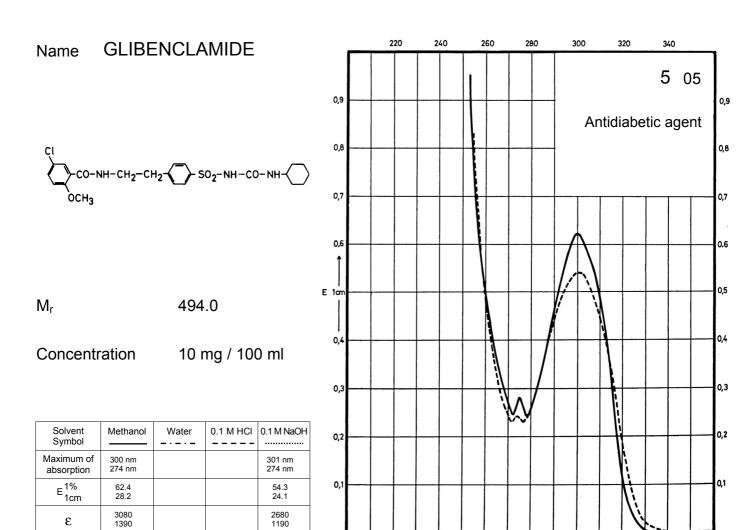


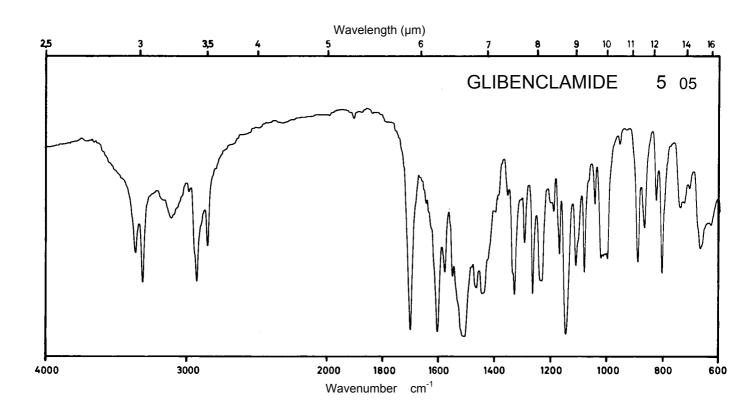

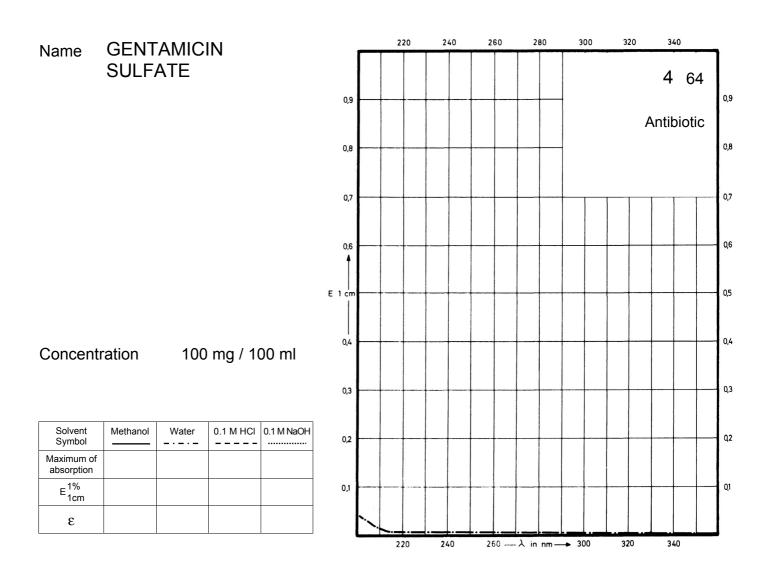


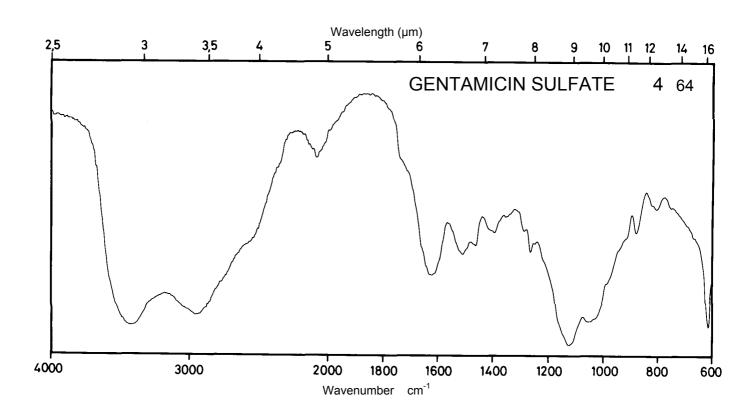

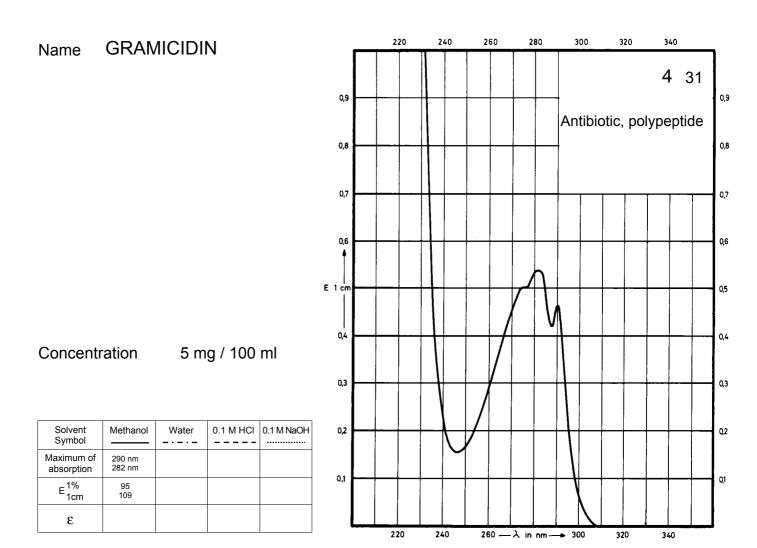

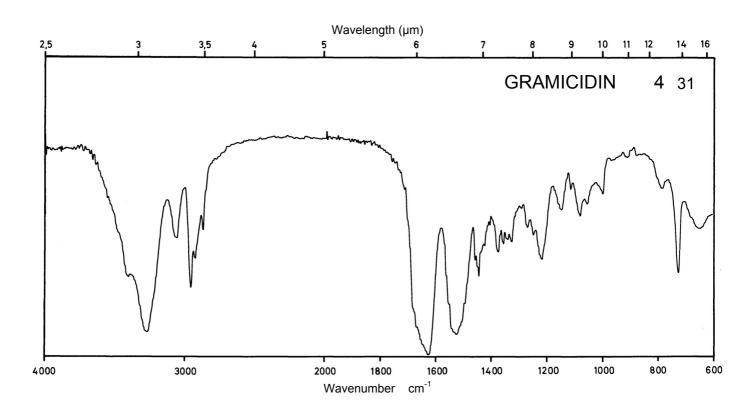





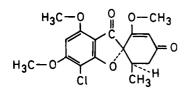

352.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	290 nm			
E ^{1%} 1cm	690			
3	24340			

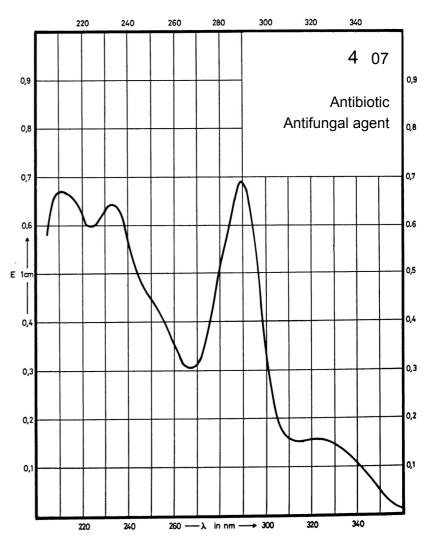

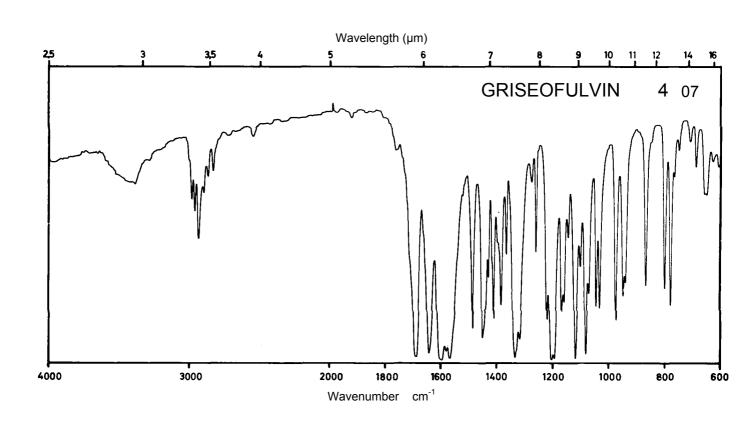



Name GLAFENINE


 M_{r}


372.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	358 nm 256 nm 228 nm		345 nm 225 nm	
E ^{1%} 1cm	522 467 1163		467 998	
ε	19460 17410 43360		17410 37210	

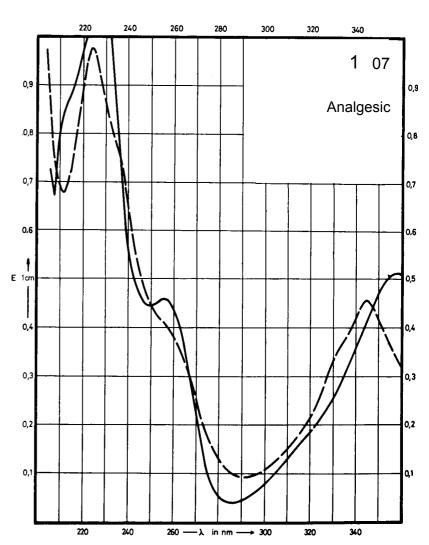


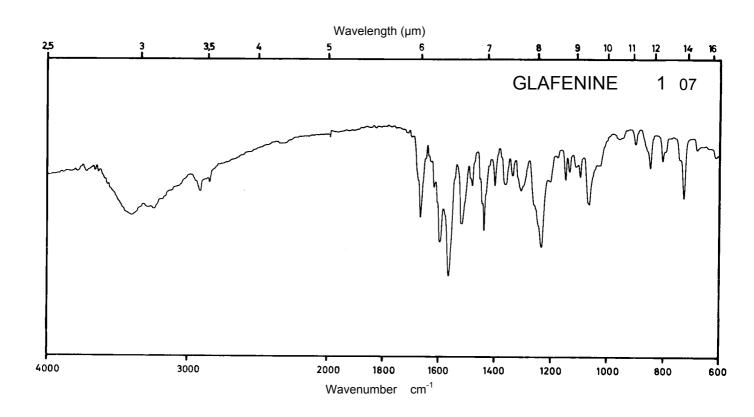

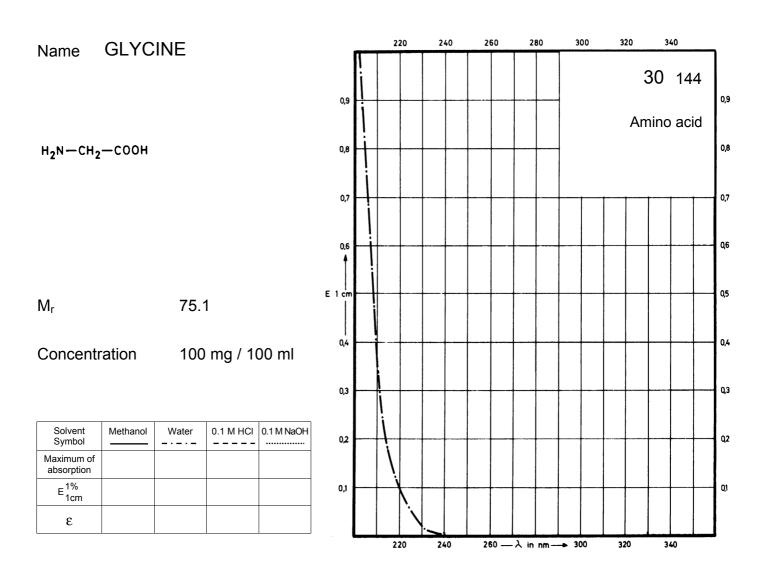

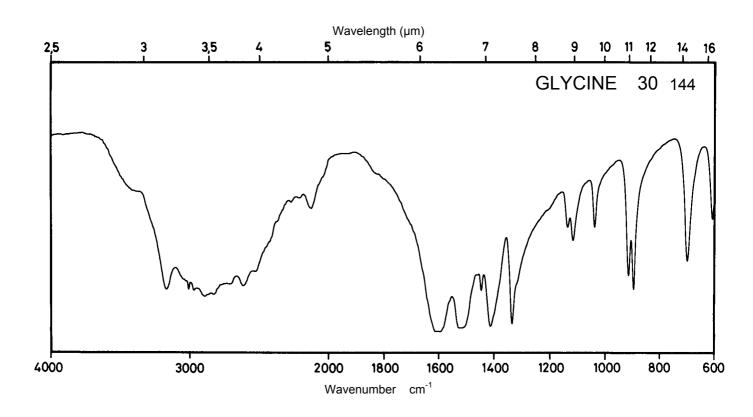





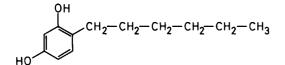


 M_r 194.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		278 nm	297 nm
E ^{1%} 1cm	151		133	226
ε	2940		2600	4400

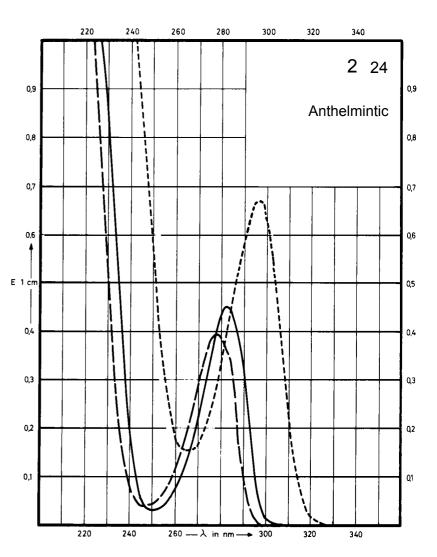

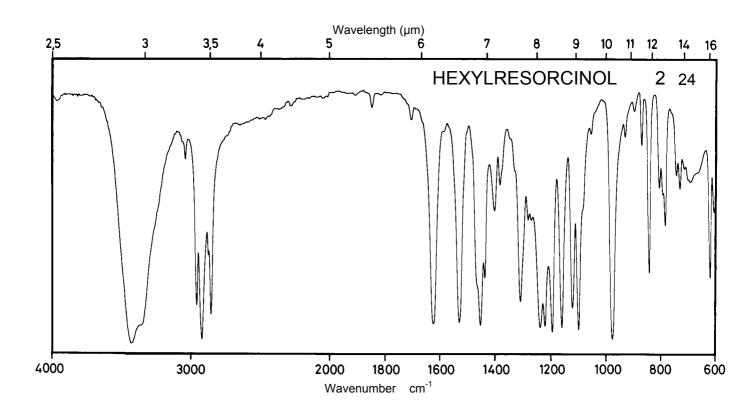


Name HISTIDINE HYDROCHLORIDE

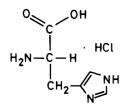


 M_{r}


191.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		211 nm	211 nm	
E ^{1%} 1cm				
3				

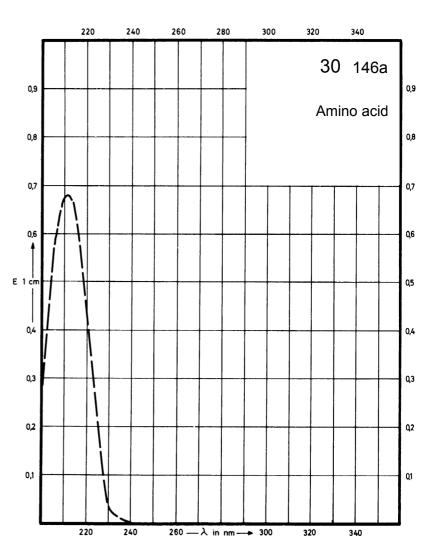

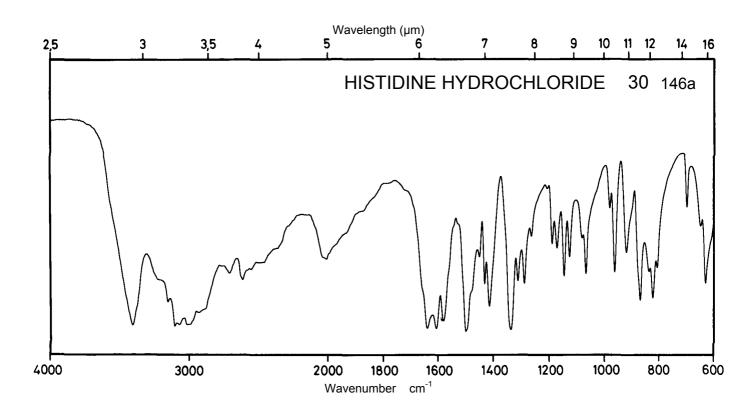




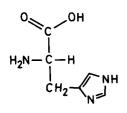


 M_r


155.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			211 nm	
E ^{1%} 1cm				
3				

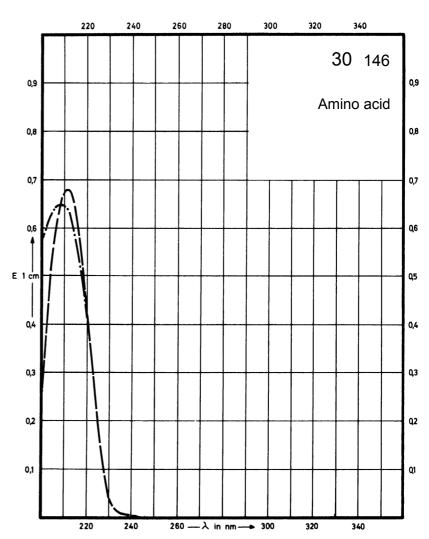

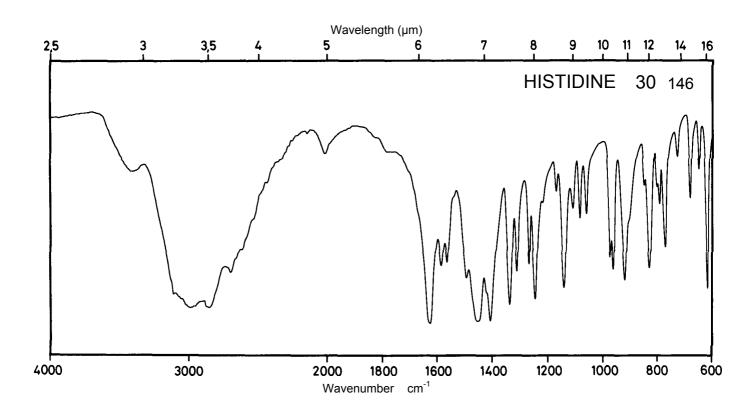




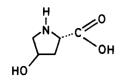


 M_r


131.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

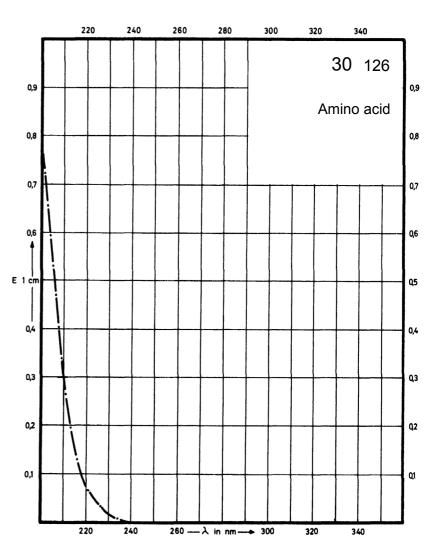

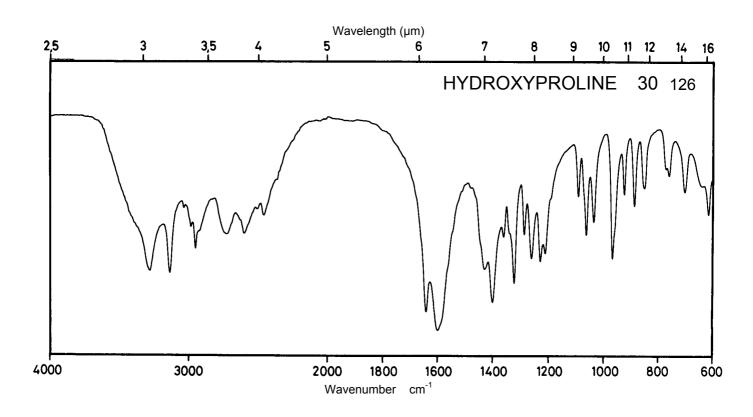

Name HYDROXYMETHYL-NICOTINAMIDE

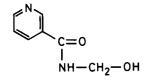



 M_r

152.2

Concentration 1.5 mg / 100 ml 3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm		261 nm	261 nm
E ^{1%} 1cm	212		377	196
ε	3230		5740	2980

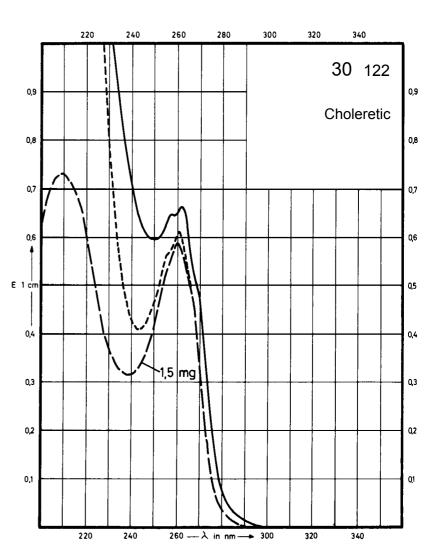

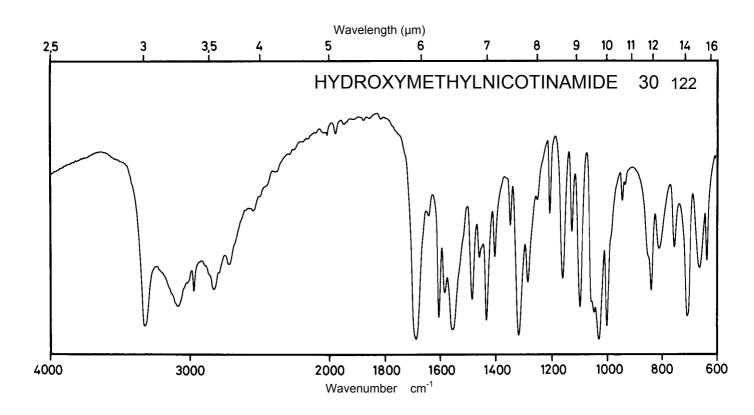



Name HALOFANTRINE HYDROCHLORIDE

 M_r


536.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm 258 nm 231 nm		310 nm 258 nm 231 nm	
E ^{1%} 1cm	261 990 526		220 865 475	
3	14000 53100 28200		11800 46400 25500	

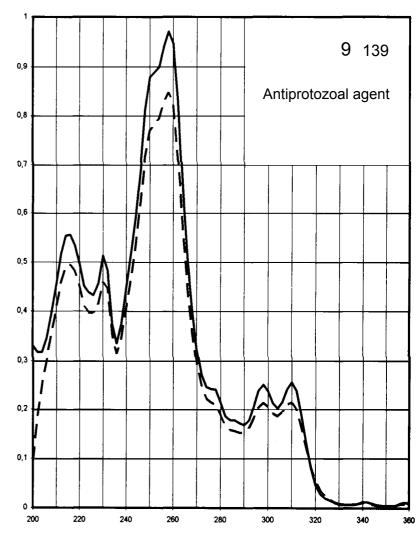

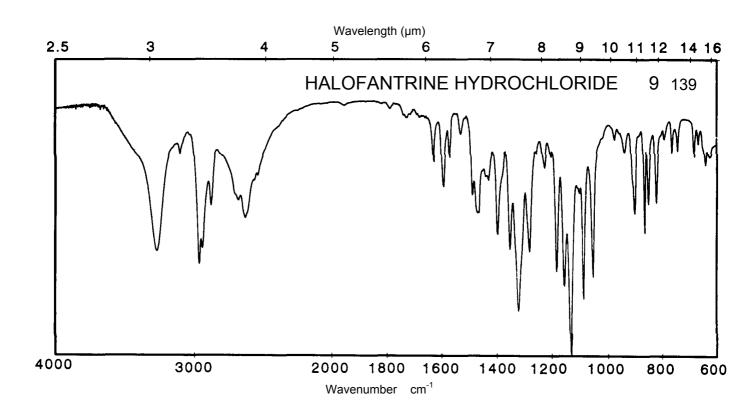




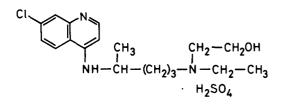


 M_{r}


434.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	344 nm 331 nm 220 nm		343 nm 256 nm 220 nm	331 nm 254 nm 219 nm
E ^{1%} 1cm	337 395 828		457 385 767	272 409 889
ε	14600 17150 35950		19800 16700 33300	11800 17750 38600

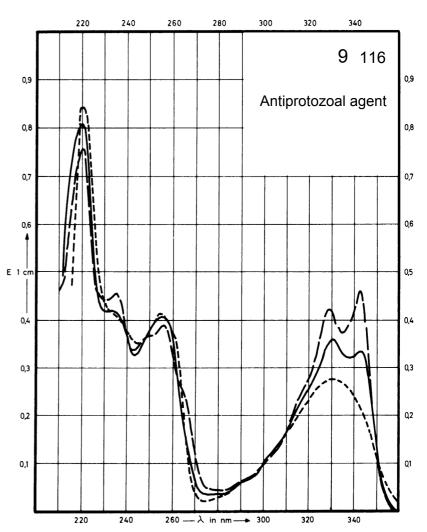

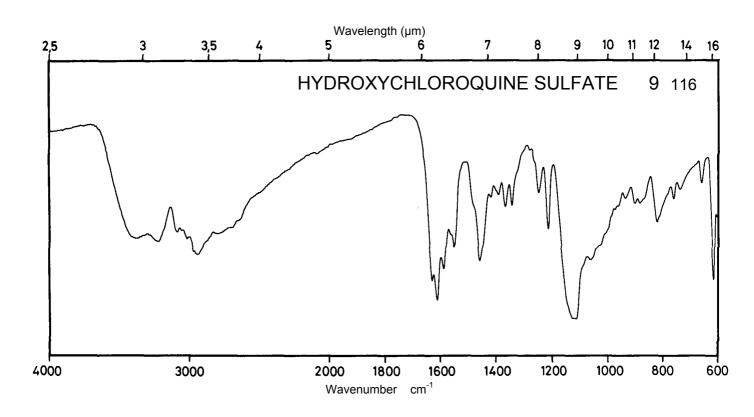





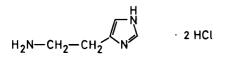


 M_r 184.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

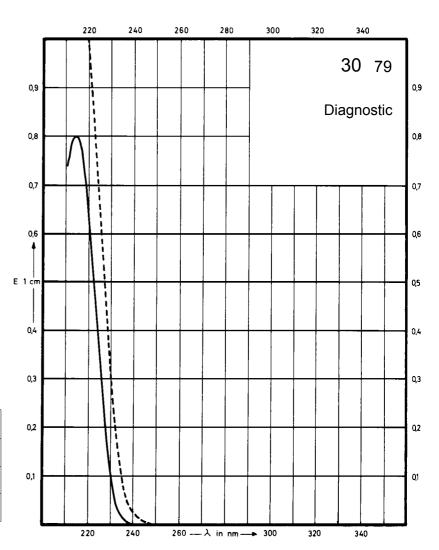

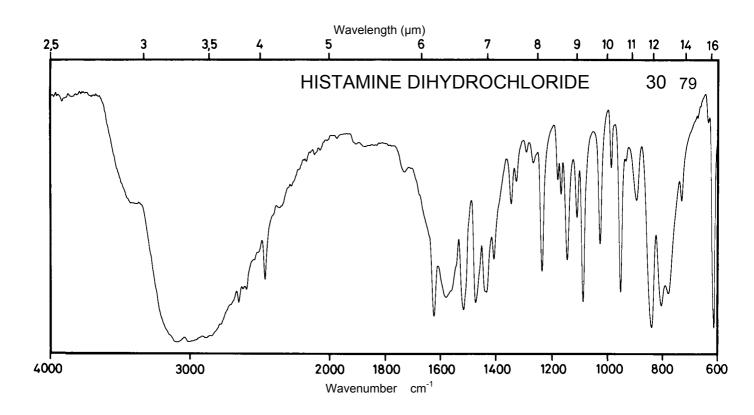



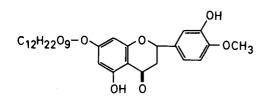


 M_{r}


610.6

Concentration


1.25 mg / 100 ml 2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	325 nm 284 nm	330 nm 284 nm	330 nm 284 nm	Decom- position observed
E ^{1%} 1cm	63 325	60 298	60 298	
ε	3850 19840	3660 18200	3660 18200	

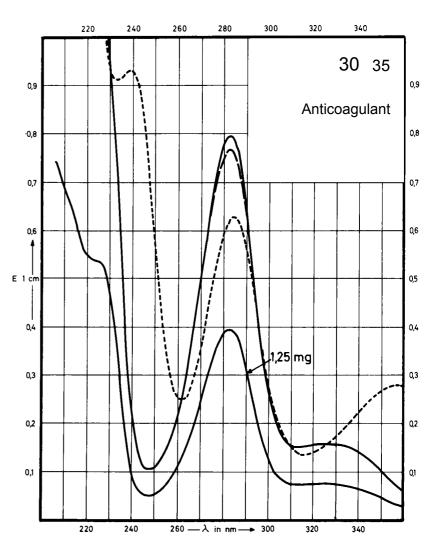

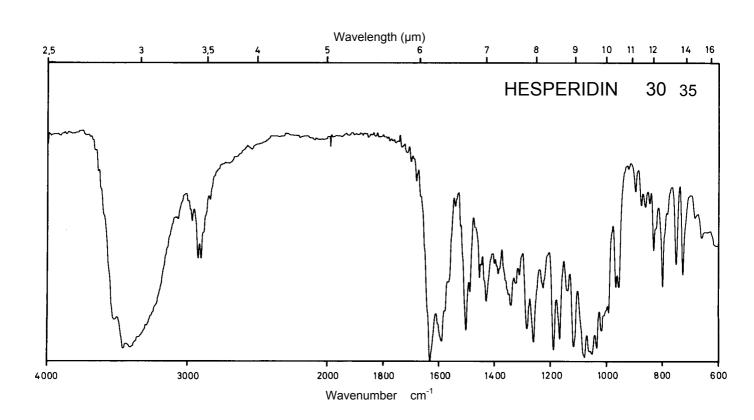


Name 4-HYDROXYPHENA-ZONE

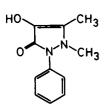


 M_{r}


204.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 252 nm	270 nm 248 nm	249 nm	Decom- position observed
E ^{1%} _{1cm}	487 578	468 479	460	
ε	9940 12000	9560 9770	9400	

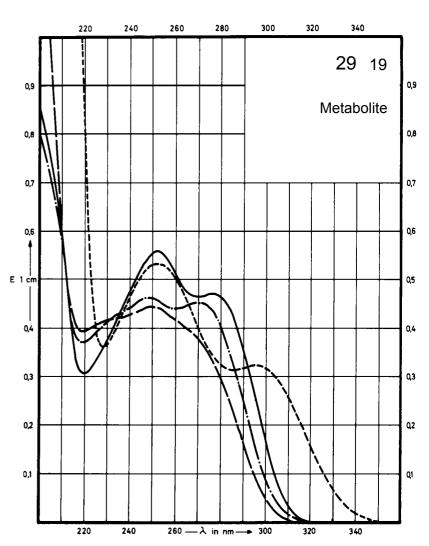

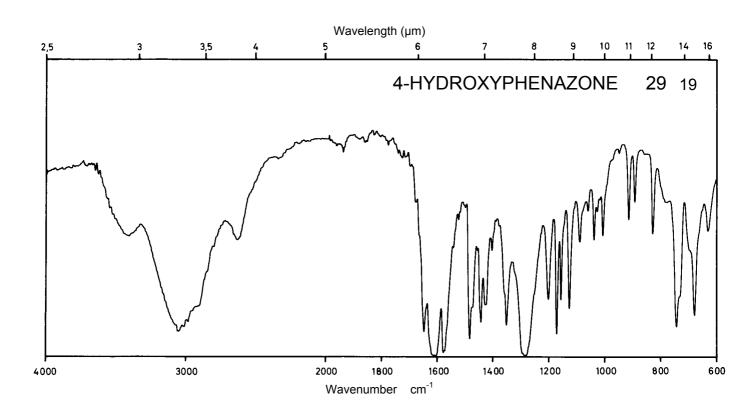



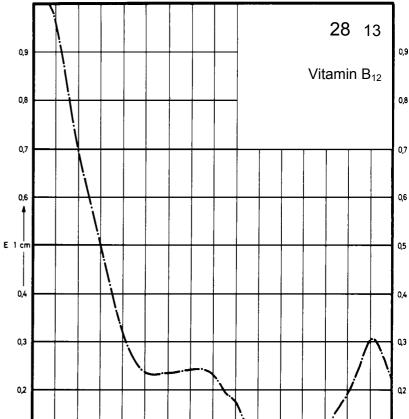

260 — λ in nm → 300

260

220

280

300

320


340

340

 M_r 1346.4


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		525 nm 351 nm		
E ^{1%} 1cm		52 155		
ε		6980 20840		

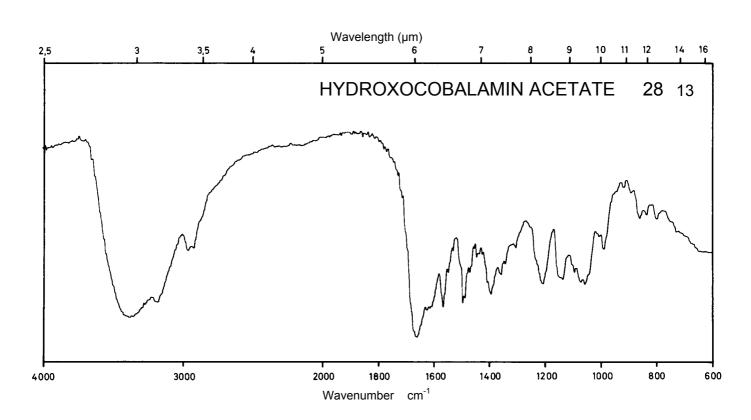

0,1


Name HOMOFENAZINE DIHYDROCHLORIDE

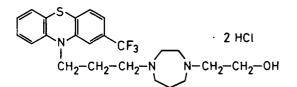


 M_{r}


524.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	309 nm 259 nm		305 nm 255 nm	
E ^{1%} 1cm	70 621		68 583	
3	3670 32600		3580 30600	

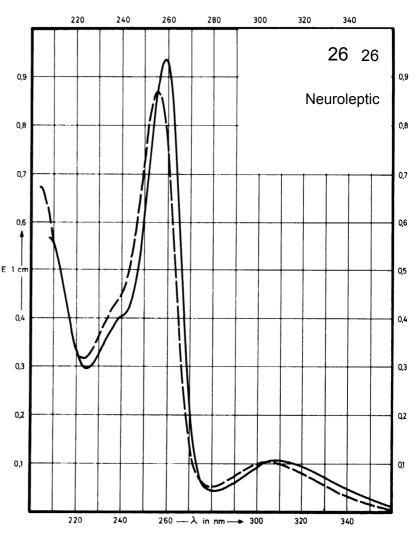

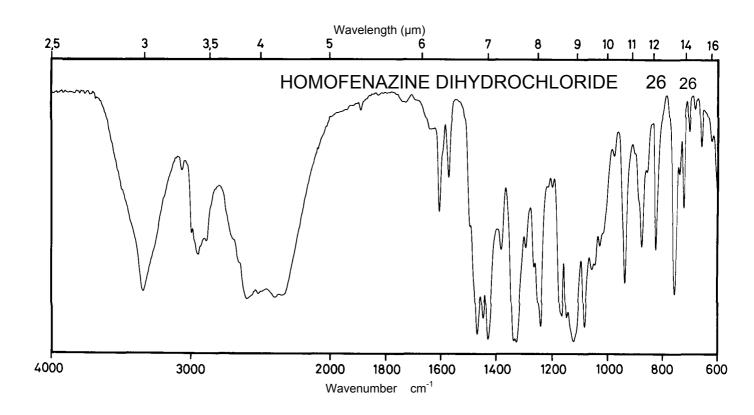



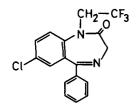


 M_{r}


352.7

Concentration

ca. 0.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	224 nm			
E ^{1%} 1cm	ca. 960			
ε	ca. 34000			

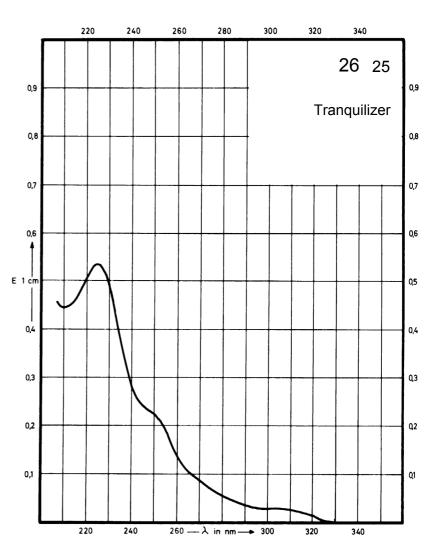


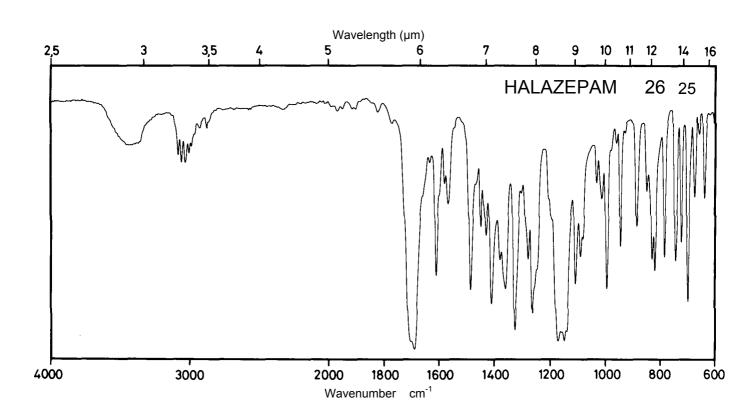


Name HYDROXYZINE DIHYDROCHLORIDE


 M_r 447.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm		270 nm 263 nm	
E ^{1%} 1cm	12.8		14.8 18.4	
ε	572		664 823	

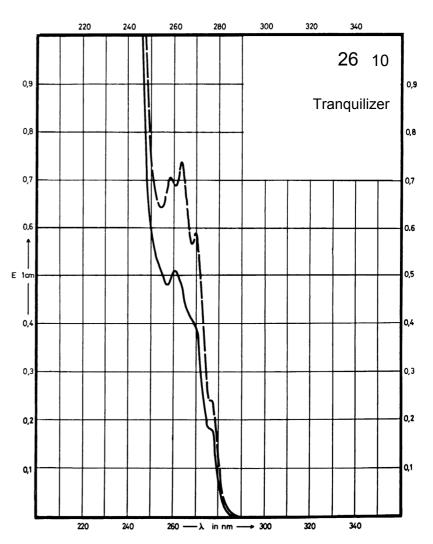

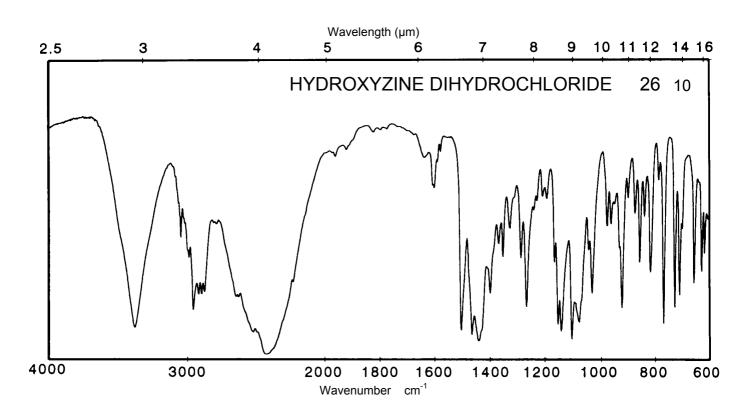





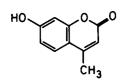

 M_r 176.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	322 nm		321 nm	Decom- position observed
E ^{1%} 1cm	867		853	
ε	15300		15000	

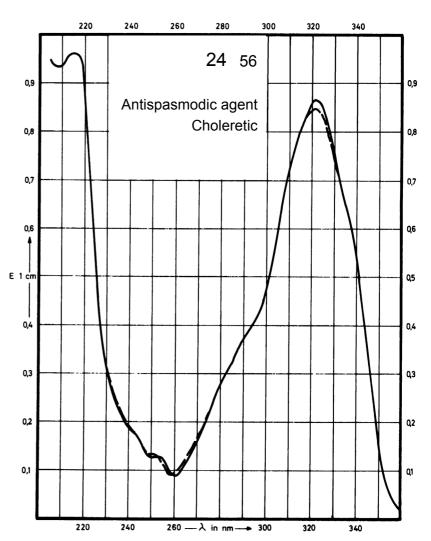

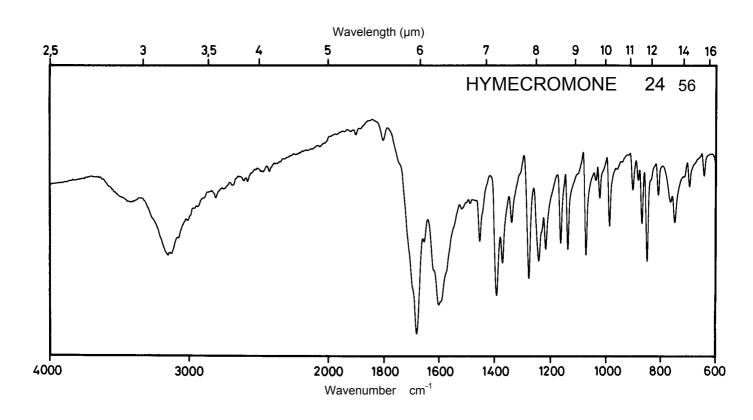


Name HEXOPRENALINE SULFATE


$$\begin{array}{c} \text{HO} \qquad \text{OH} \\ \text{HO} \qquad \text{CH-CH}_2\text{-NH} \\ \text{HO} \qquad \text{CH-CH}_2\text{-NH} \\ \text{OH} \end{array}$$

 M_{r}


518.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 224 nm	280 nm 221 nm	280 nm 221 nm	Decom- position observed
E ^{1%} 1cm	118 222	105 230	106 233	
3	6100 11500	5400 11900	5500 12000	

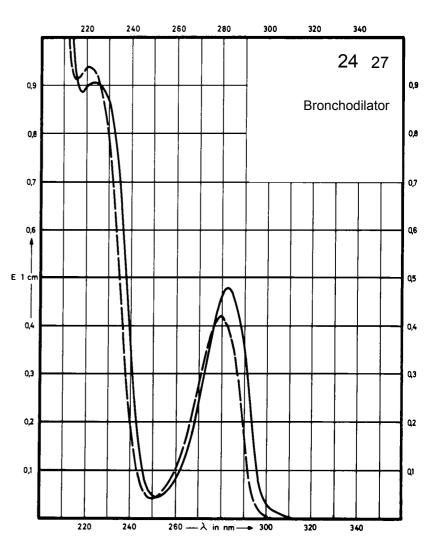

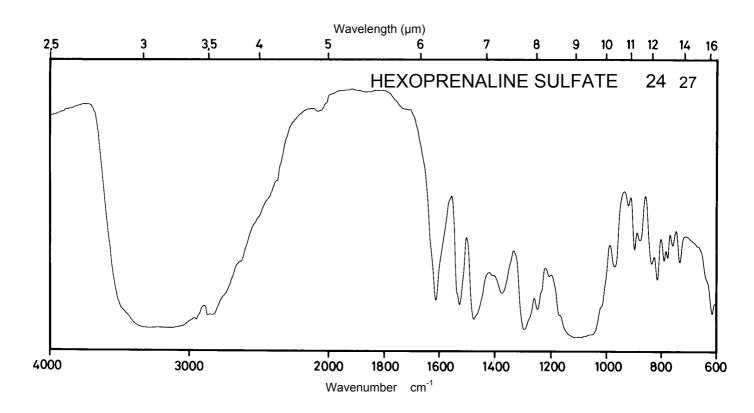




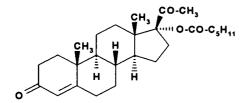


 M_r


428.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm			
E ^{1%} 1cm	402			
ε	17200			

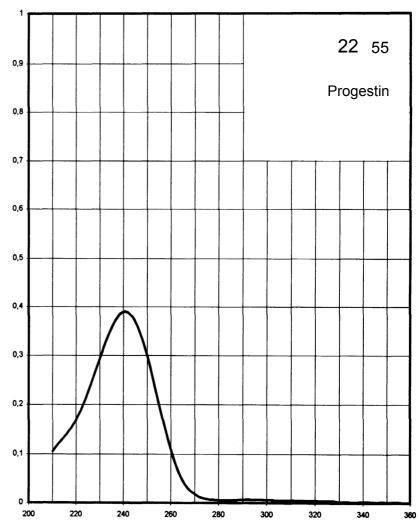

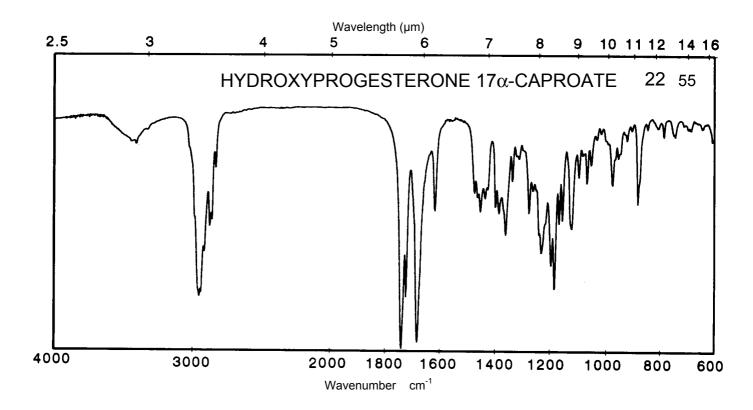

Name HYDROXYPRO-GESTERONE

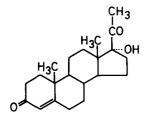


 M_{r}


330.5

Concentration


1.2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	242 nm	249 nm		
E ^{1%} 1cm	500	492		
ε	16530	16260		

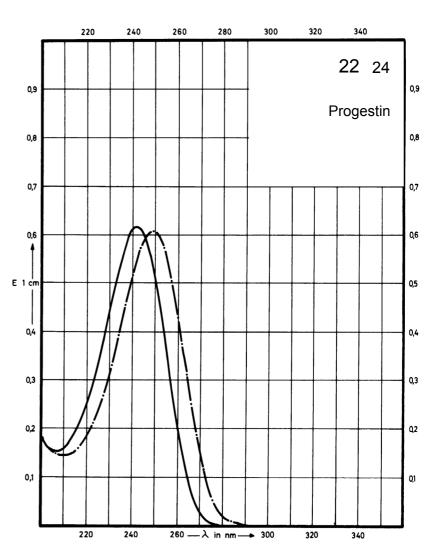

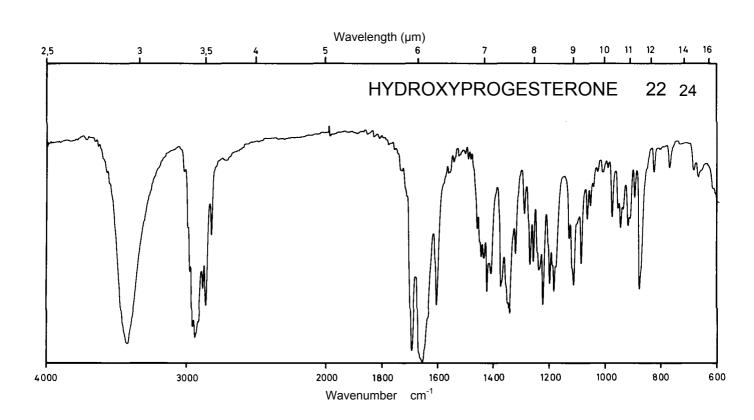




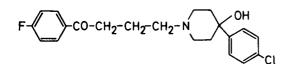


 M_{r}

375.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm		248 nm	
E ^{1%} 1cm	84		86	
ε	3160		3230	

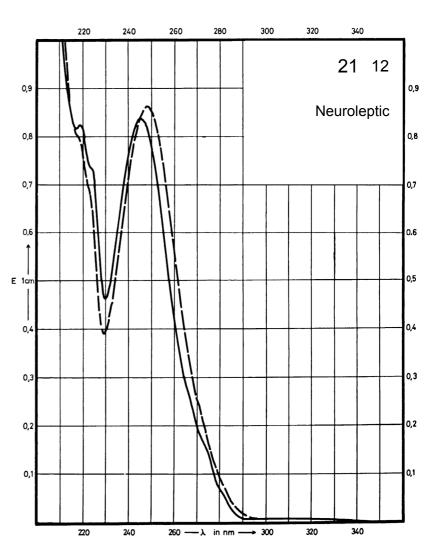

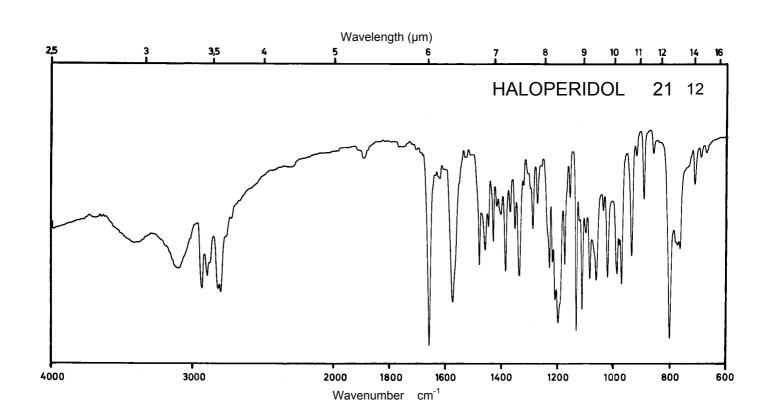





 M_{r}

356.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm		257 nm	
E ^{1%} 1cm	5.6		6.0	
ε	200		210	

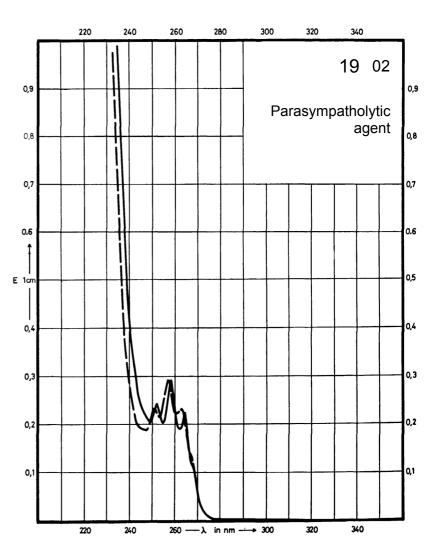

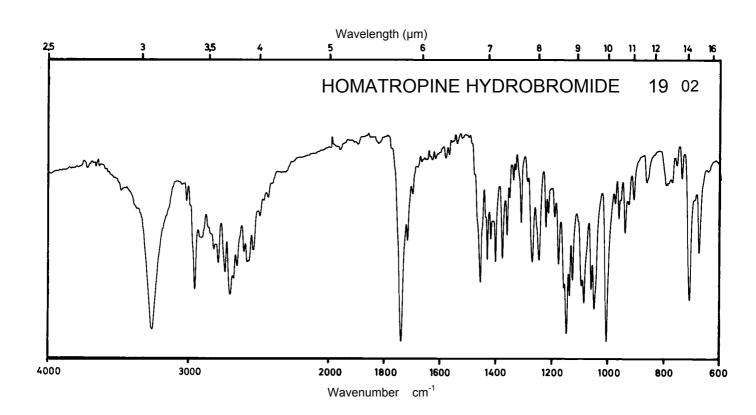


Name HYDROXYPROCAINE HYDROCHLORIDE


$$\begin{bmatrix} \mathsf{OH} & \mathsf{CH}_2 - \mathsf{CH}_3 \\ \mathsf{H}_2 \mathsf{N} - \mathsf{COO} - \mathsf{CH}_2 - \mathsf{CH}_2 - \mathsf{N} + \mathsf{CH}_2 - \mathsf{CH}_3 \\ \mathsf{CH}_2 - \mathsf{CH}_3 \end{bmatrix} \mathsf{CI}^-$$

 M_r

288.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	307 nm	305 nm	303 nm	325 nm
	242 nm	235 nm	236 nm	273 nm
E ^{1%}	740	605	350	288
1cm	295	227	300	377
ε	21370	17470	10110	8320
	8520	6560	8560	10890

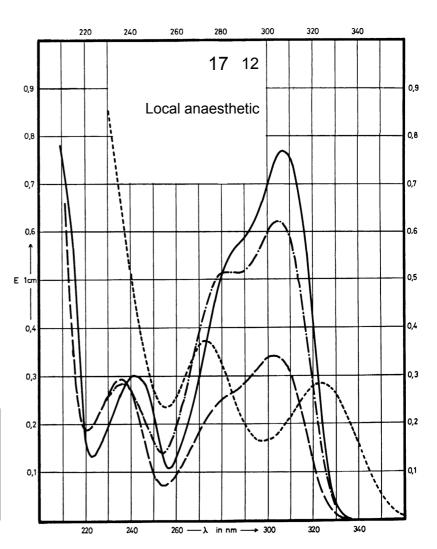

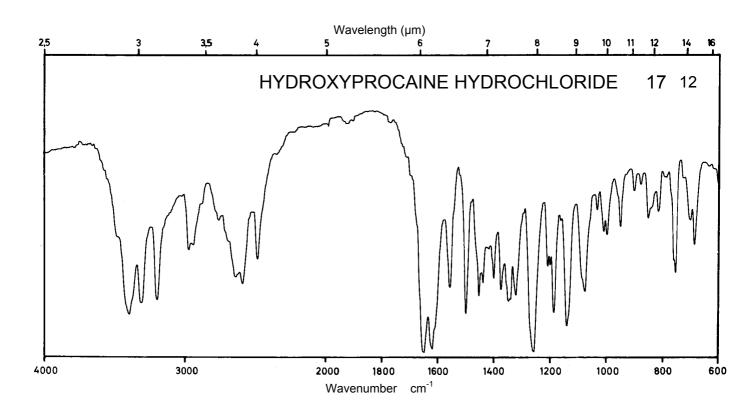



Name HYDROXYTETRACAINE HYDROCHLORIDE

 M_{r}


316.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm 248 nm		310 nm 237 nm	322 nm 285 nm
E ^{1%} 1cm	980 257		286 294	410 405
ε	31060 8140		9060 9320	12990 12830

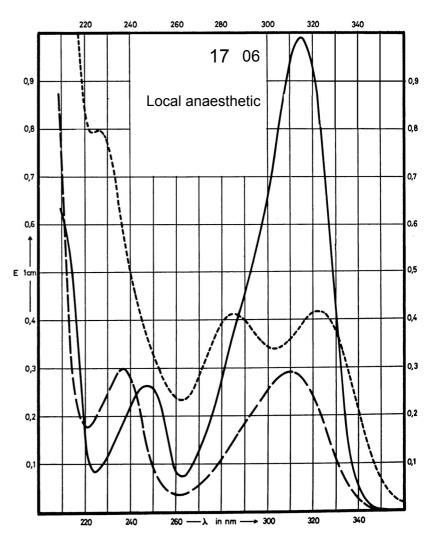

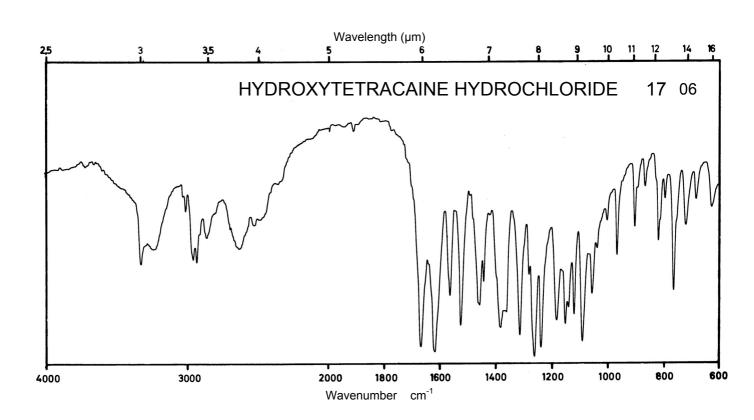





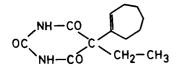


 M_r 250.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				251 nm
E ^{1%} 1cm				274
3				6860

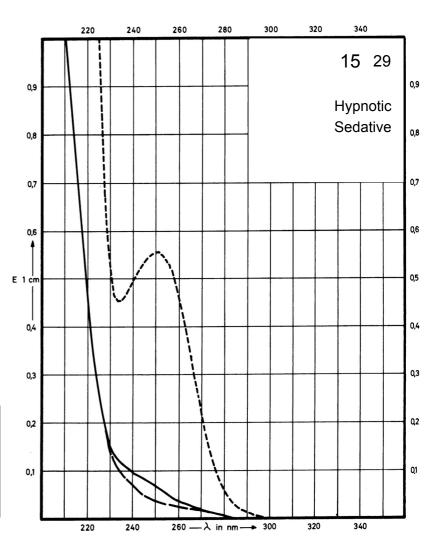

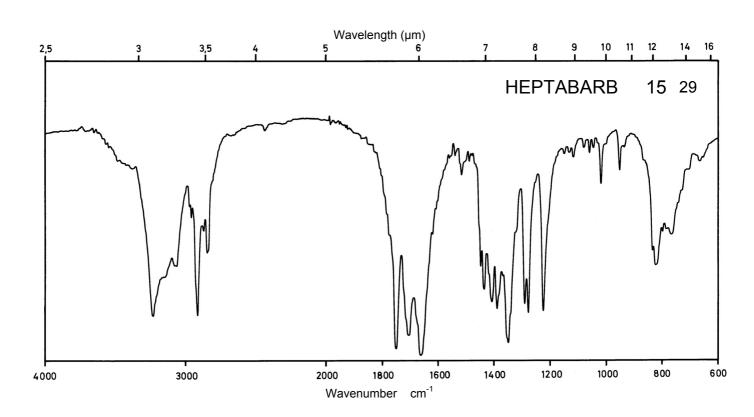




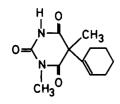



 M_{r}


236.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				243 nm
E ^{1%} 1cm				328
ε				7750

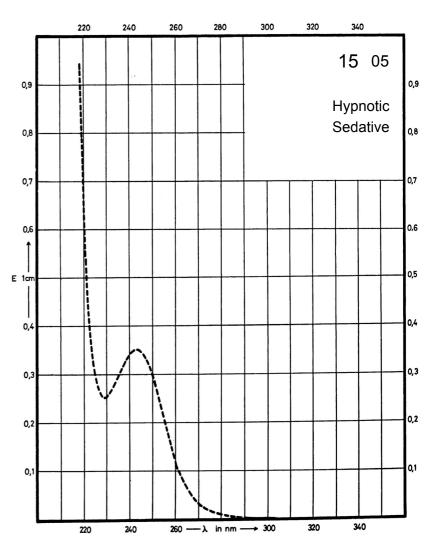

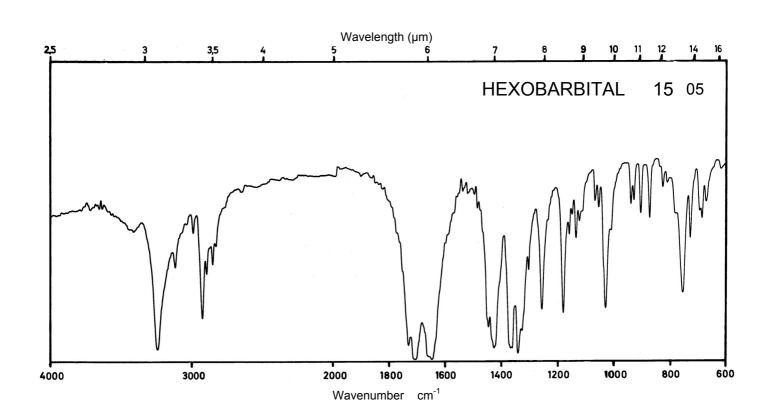

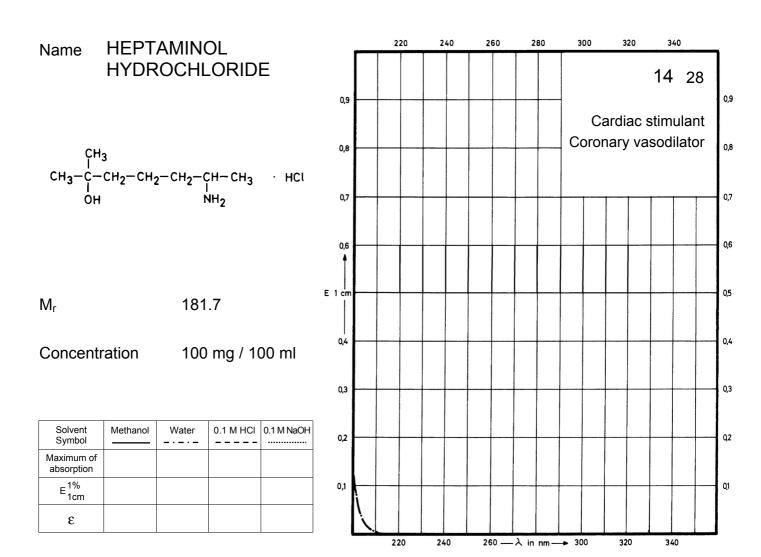

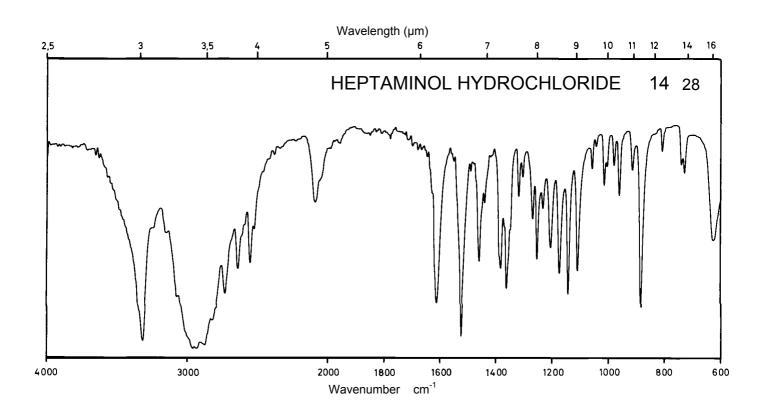



Name HEXOBENDINE DIHYDROCHLORIDE

 M_{r}

665.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm		266 nm	Decom- position observed
E ^{1%} 1cm	303		286	
3	20190		19010	

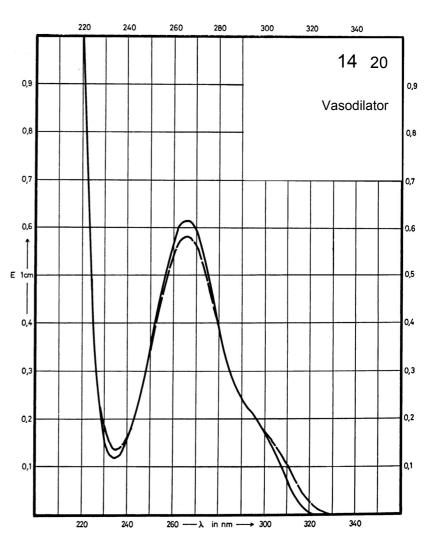

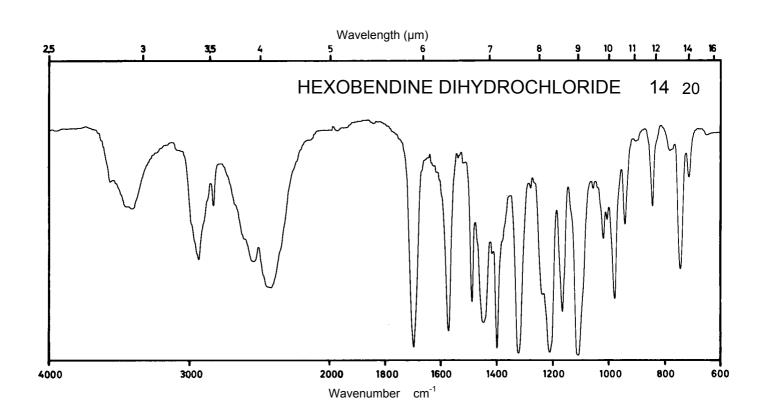



Name HYDROFLUME-THIAZIDE

 M_{r}


331.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 324 nm		273 nm 324 nm	273 nm 330 nm
E ^{1%} 1cm	600 118		576 122	470 102
3	19880 3910		19080 4040	15570 3380

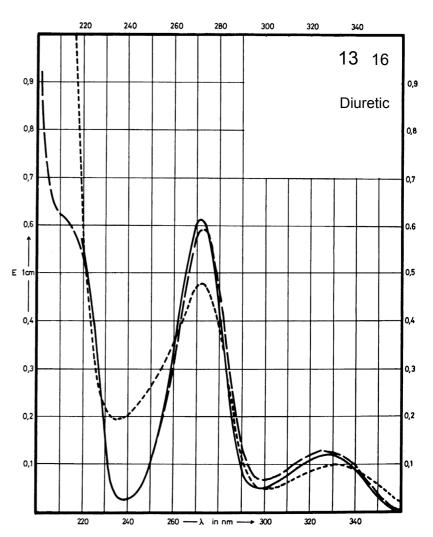

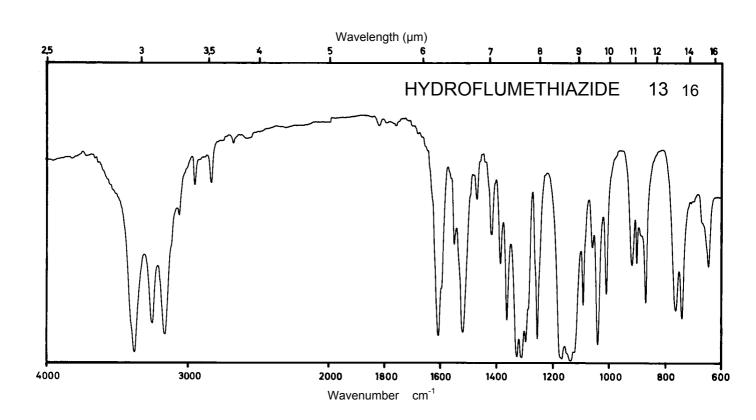


Name HYDROCHLORO-THIAZIDE

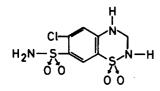


 M_{r}


297.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm 269 nm		315 nm 271 nm	320 nm 272 nm
E ^{1%} 1cm	99 668		99 651	91 523
3	2950 19890		2950 19390	2710 15570

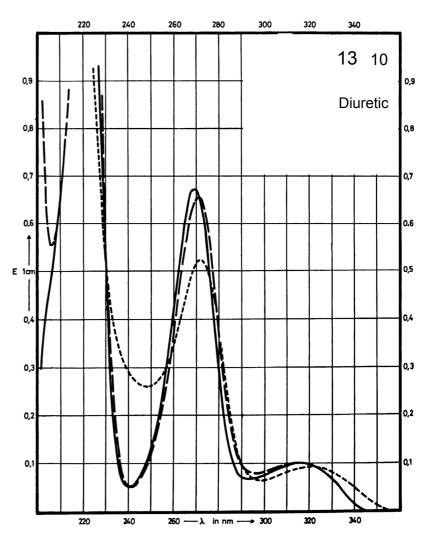

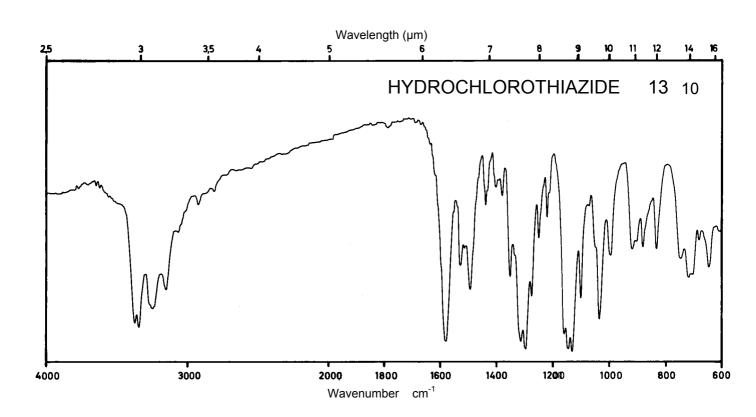





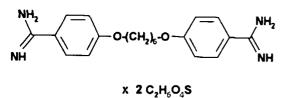


 M_{r} 606.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm	262 nm	262 nm	
E ^{1%} 1cm	551	489	492	
3	33400	29700	29900	

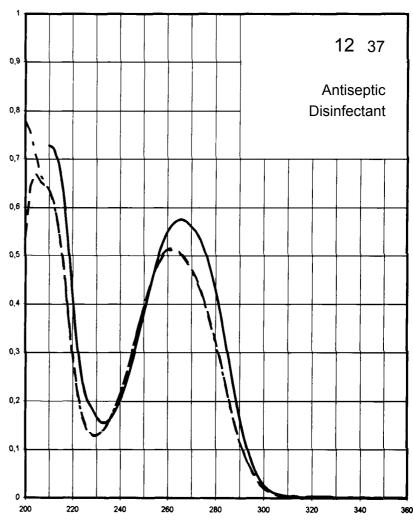

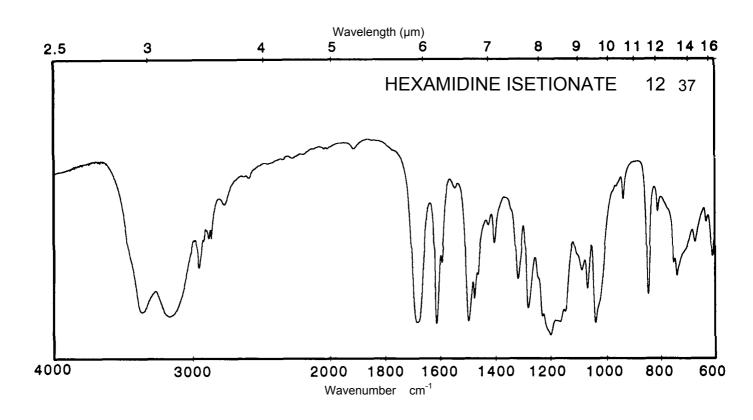





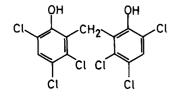


 M_{r} 406.9

Concentration 1 mg / 100 ml 5 mg / 100 ml

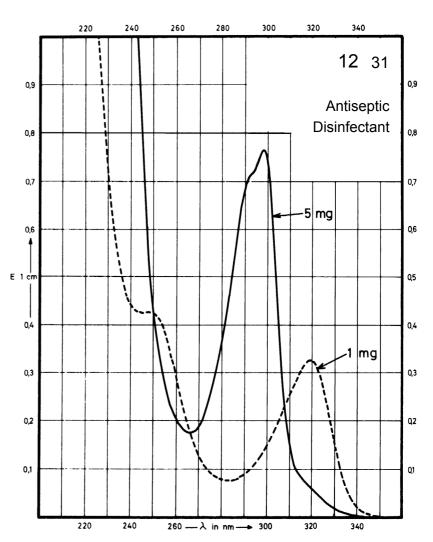

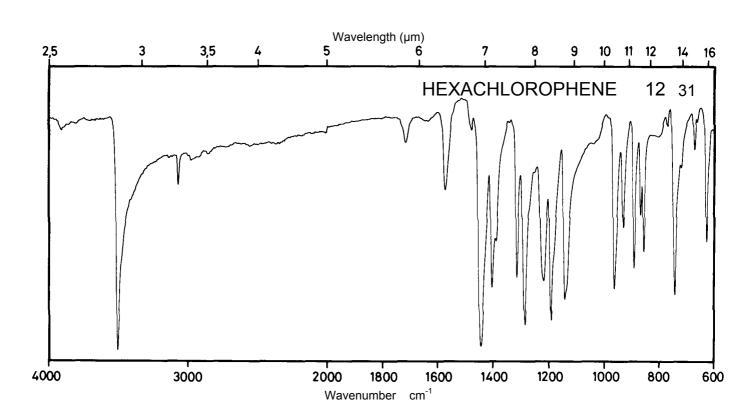




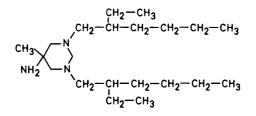


 M_{r} 339.6

Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

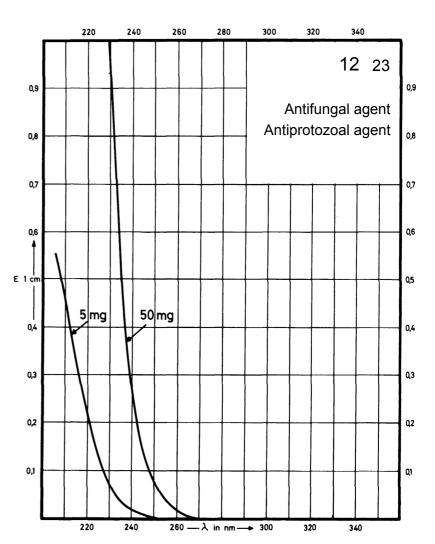

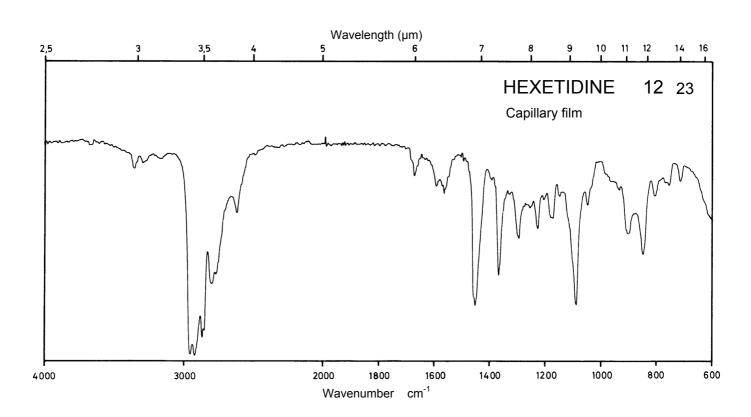

Name HALQUINOL

CI N


 M_{r}


214.1

Concentration

0.4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	330 nm 249 nm		366 nm 259 nm	343 nm 259 nm
E ^{1%} _{1cm}	141 1644		120 2234	198 1480
ε	3020 35200		2570 47830	4240 31690

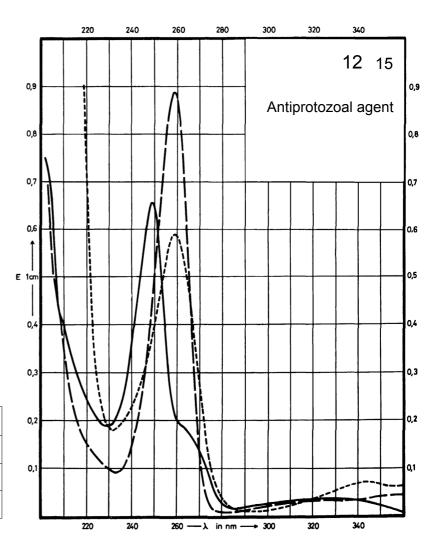

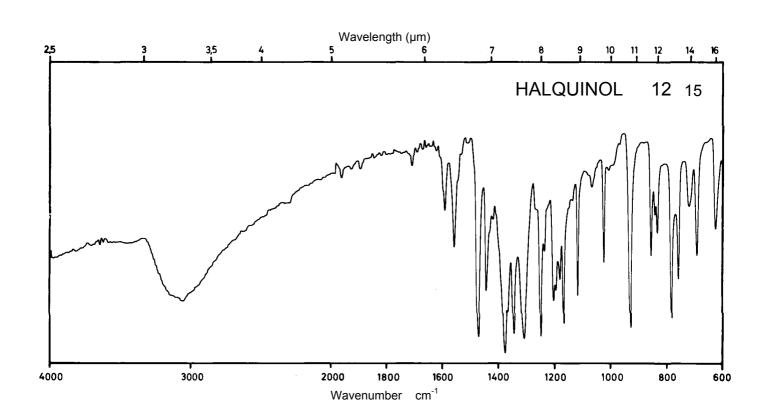

Name 4-HYDROXYBENZOIC ACID




 M_{r}


138.1

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	252 nm		255 nm	279 nm
E ^{1%} 1cm	1068		1078	1290
ε	14750		14890	17810

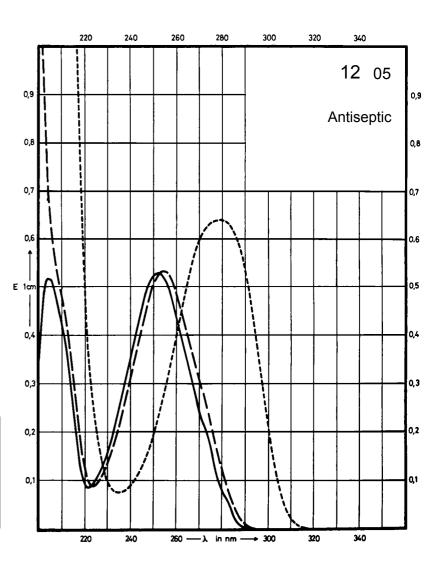

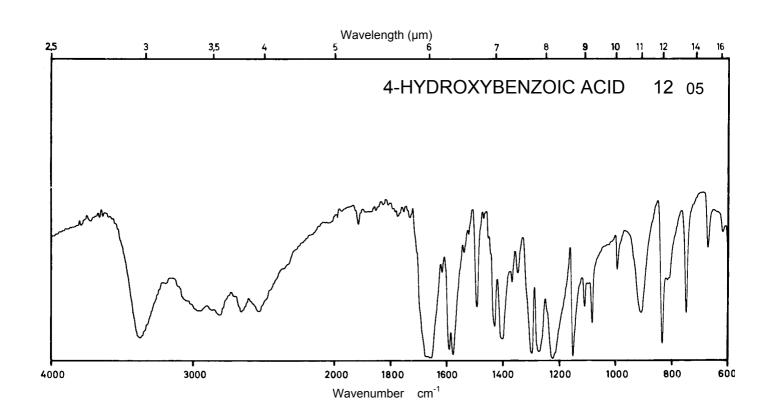





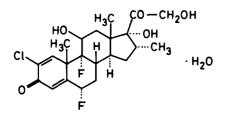


 M_{r} 462.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm	249 nm		
E ^{1%} 1cm	340	345		
ε	15730	16000		

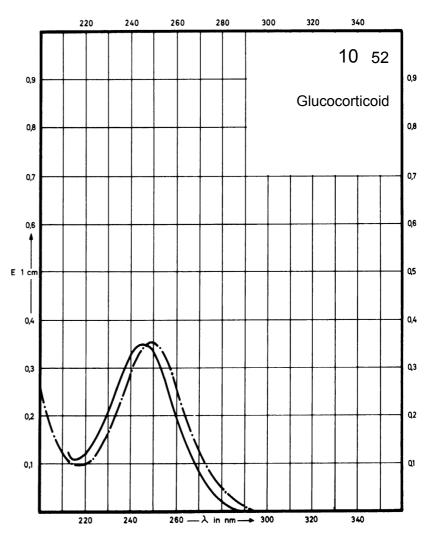

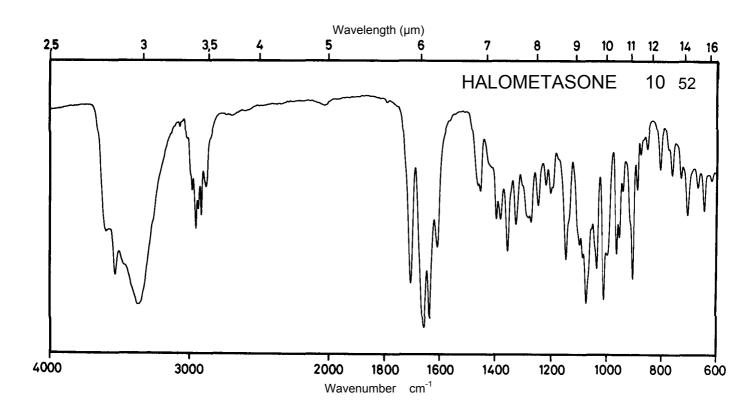



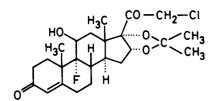



455.0

Concentration

1.05 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm			
E ^{1%} _{1cm}	377			
ε	17170			

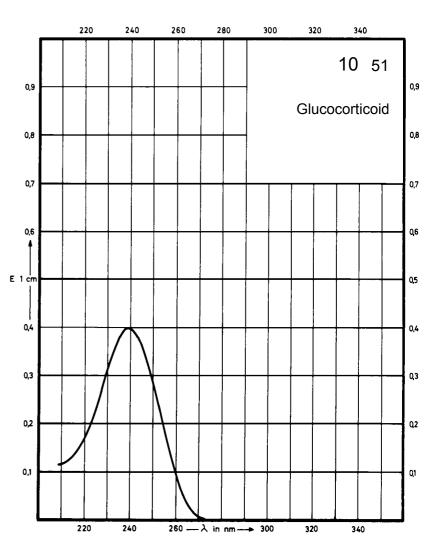

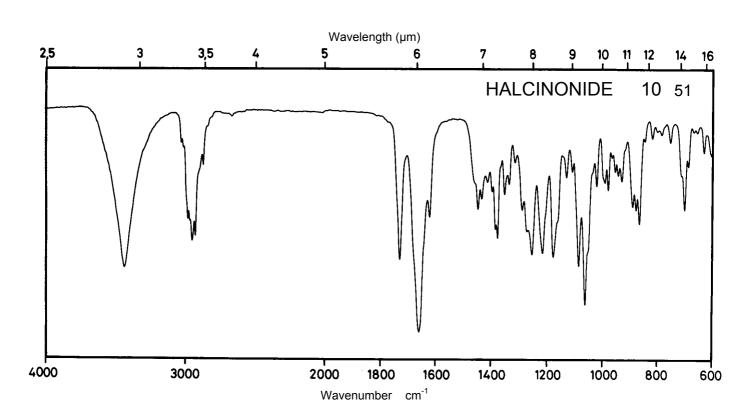


Name HYDROCORTISONE 21-ACETATE 17-PROPIONATE


 M_{r}

460.6


Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm	247 nm		
E ^{1%} 1cm	358	361		
3	16490	16650		

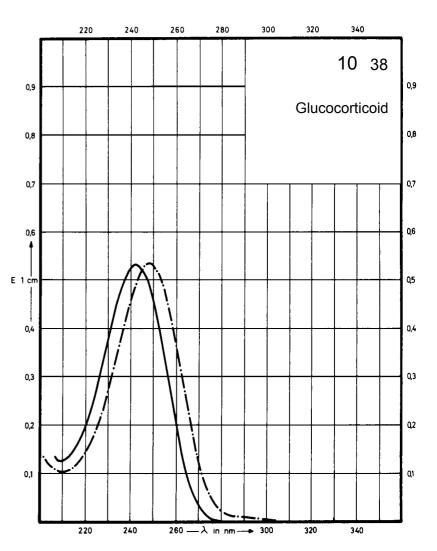

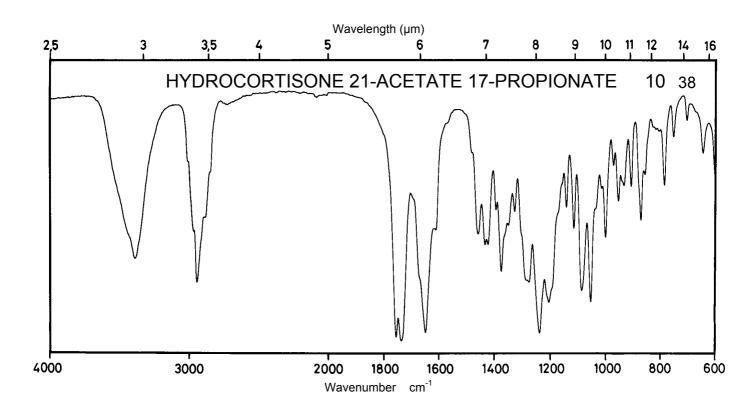




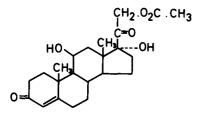


 M_r 404.5


Concentration 1.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm			
E ^{1%} 1cm	391			
ε	15820			

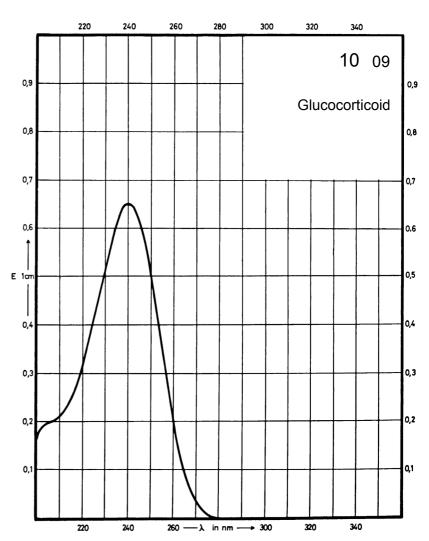

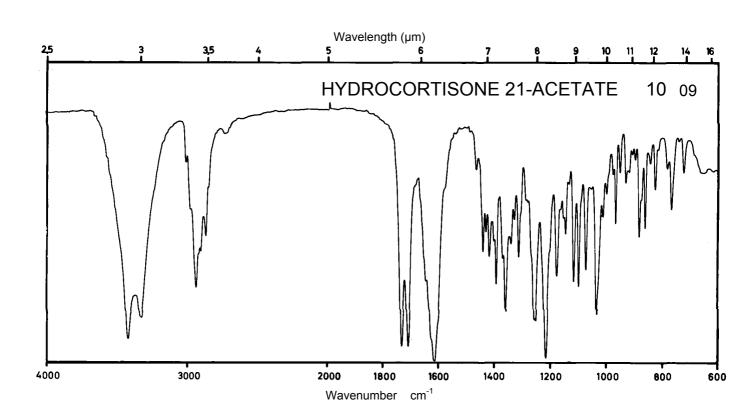



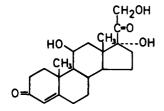



 M_r


362.5


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm			
E ^{1%} 1cm	441			
ε	15990			

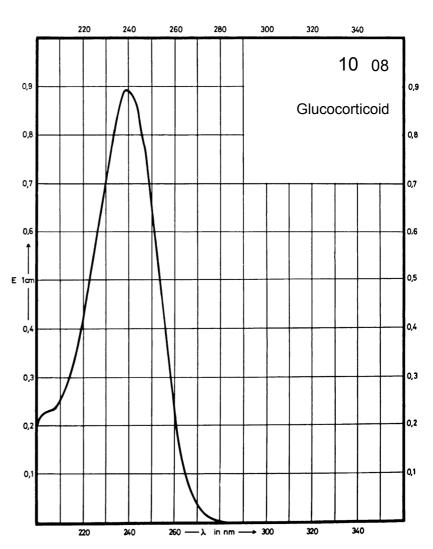


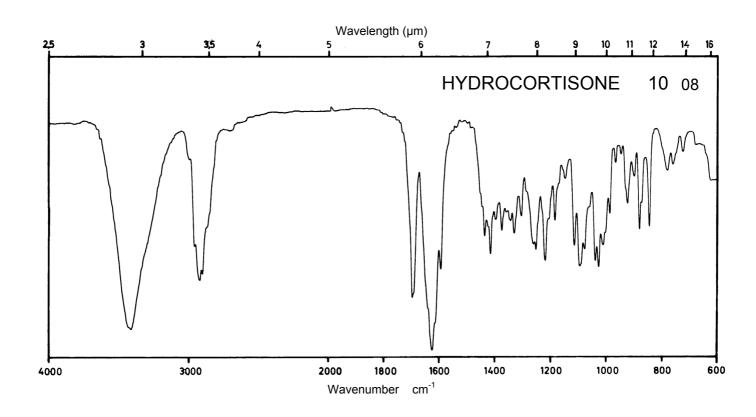

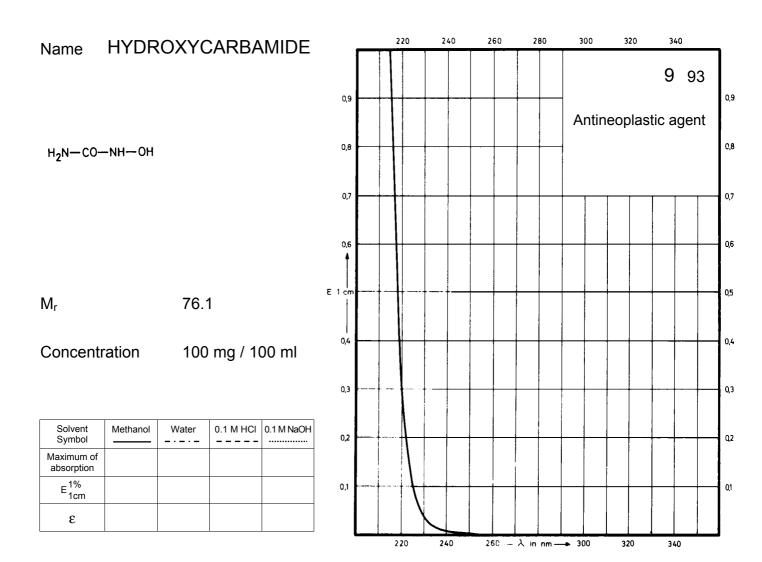

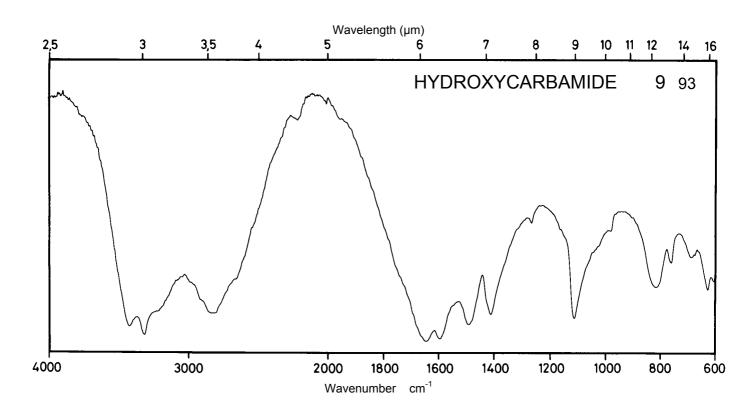




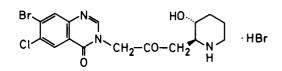


 M_{r} 495.6


Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of	326 nm 275 nm		326 nm 282 nm	
absorption	243 nm		243 nm	
_1%	54		56	
E ^{1%} 1cm	168		137	
	925 2690		1016 2780	
3	8300		6790	
C	45840		50350	

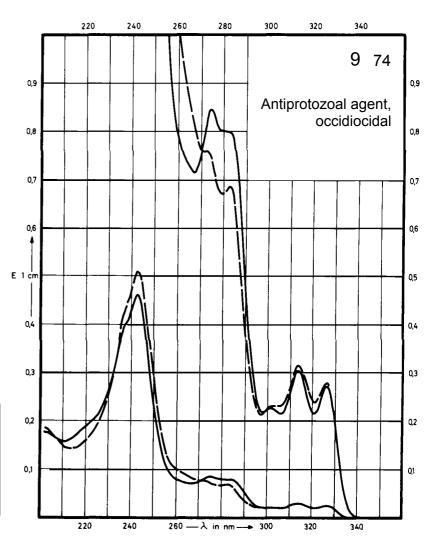

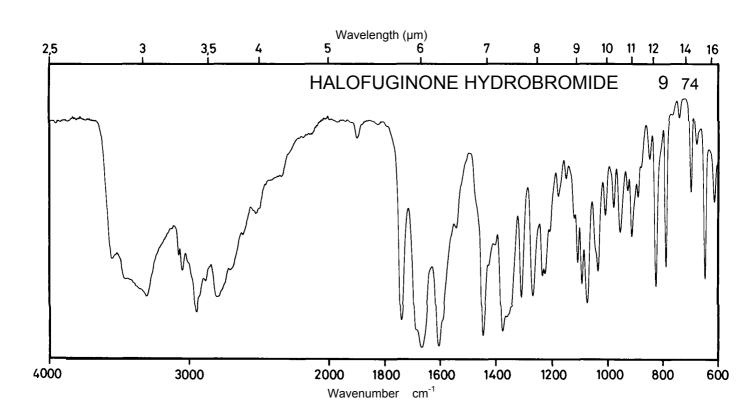


Name HYDROCODONE HYDROCHLORIDE

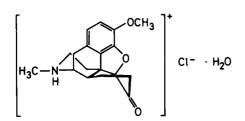


 M_r 353.8


Concentration 1.5 mg / 100 ml 15 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm		280 nm	280 nm
E ^{1%} 1cm	38.3		35.1	32.6
ε	1350		1240	1150

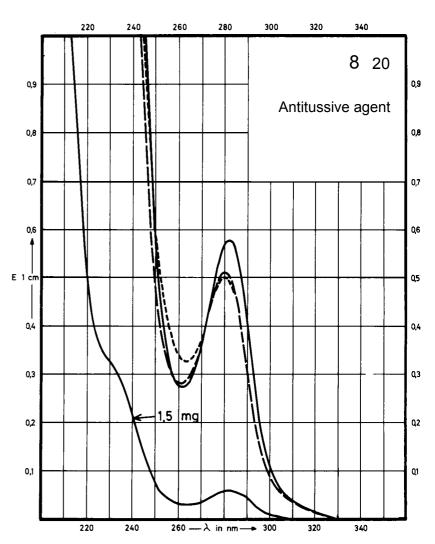

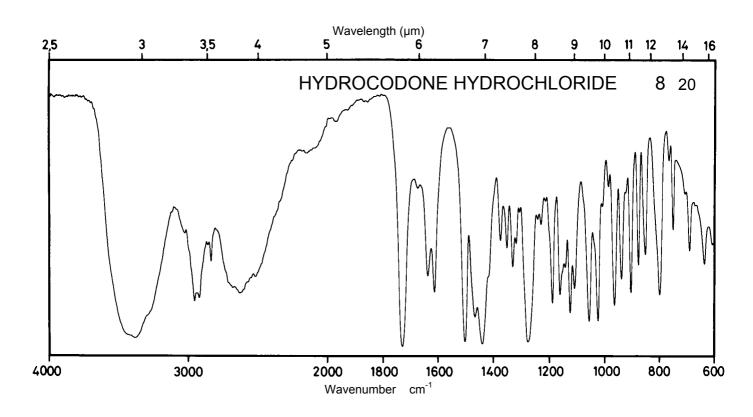




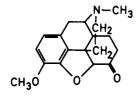


 M_{r} 299.4


Concentration 17.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm			
E ^{1%} 1cm	40			
3	1200			

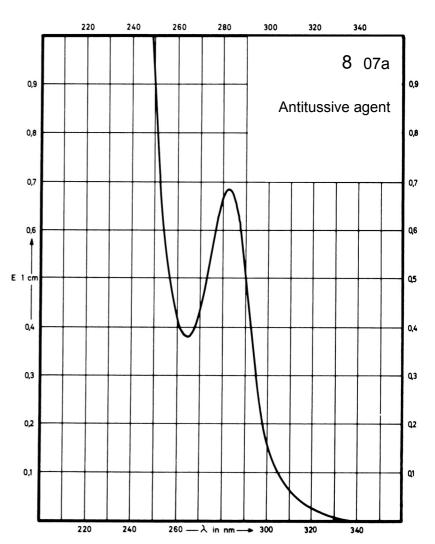

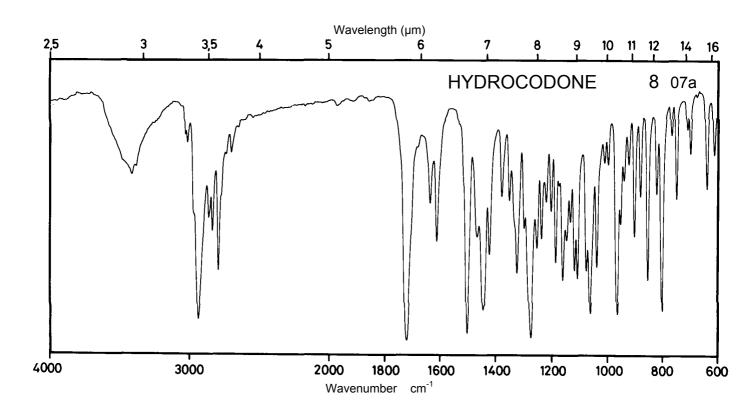




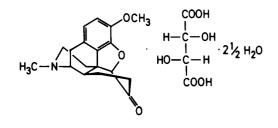


 M_{r} 494.5

Concentration 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm		280 nm	280 nm
E ^{1%} _{1cm}	28.0		26.6	24.6
ε	1380		1320	1210

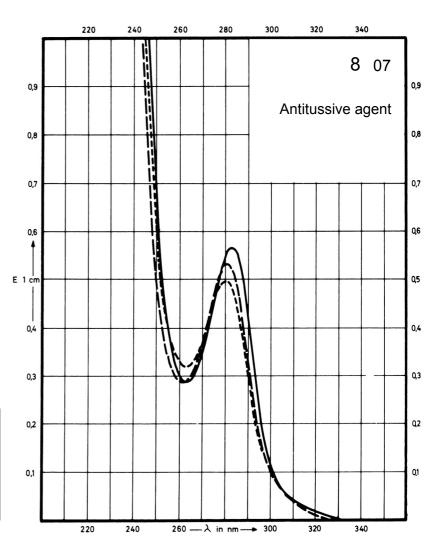

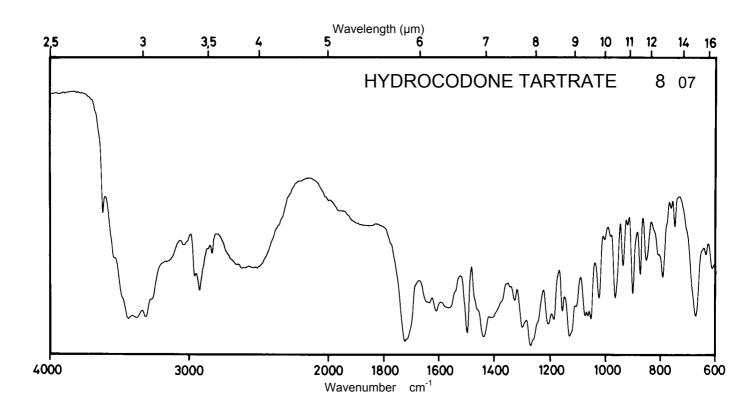



Name HYDRALAZINE SULFATE


 M_{r} 418.4


Concentration 1 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 317 nm 305 nm 262 nm 237 276 315 nm 303 nm 260 nm 233 275 Decom-position observed Maximum of absorption E^{1%} 1cm 522 9910 500 9750 11540 21850 11520 3 20940

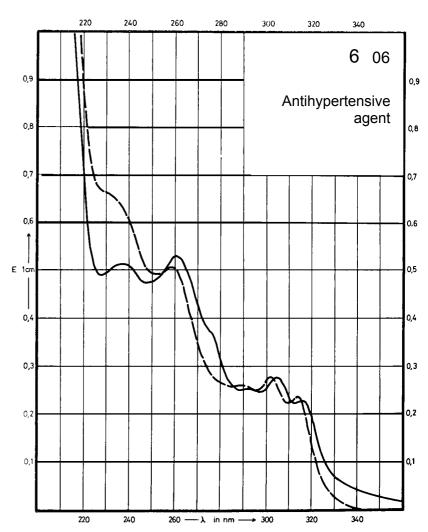


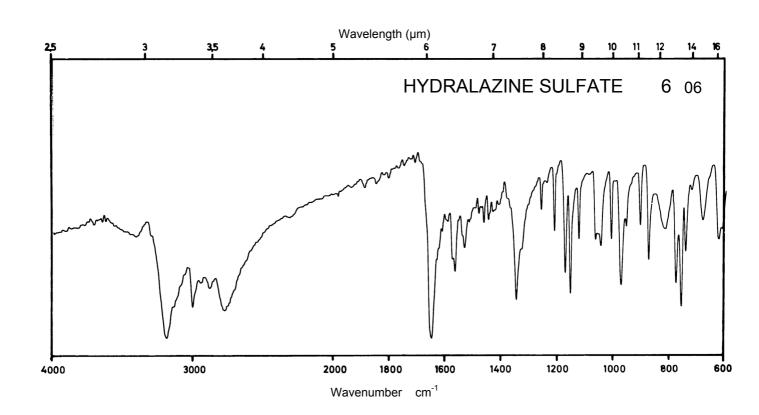

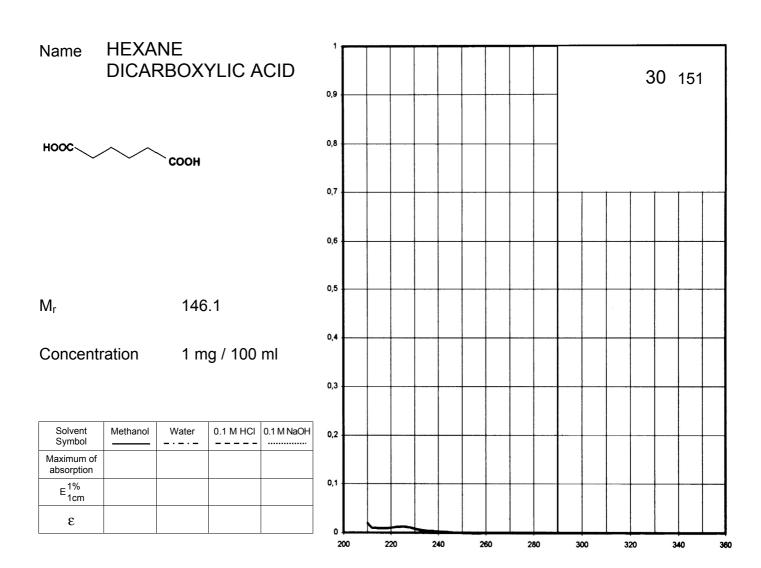

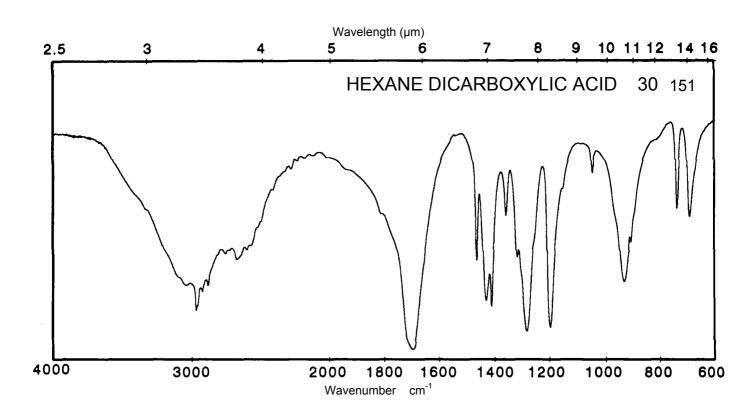



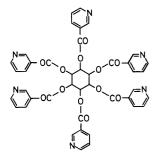


 M_r


810.7

Concentration


1.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			261 nm	262 nm
E ^{1%} 1cm			369	232
ε			32100	18800

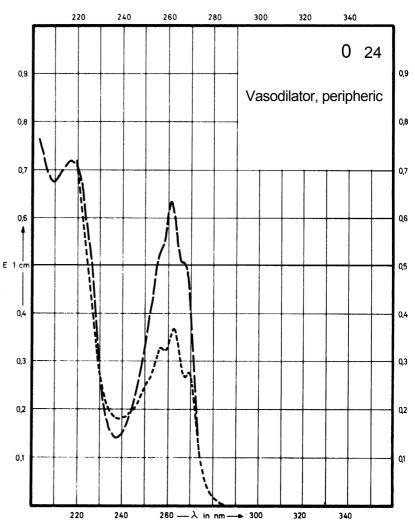

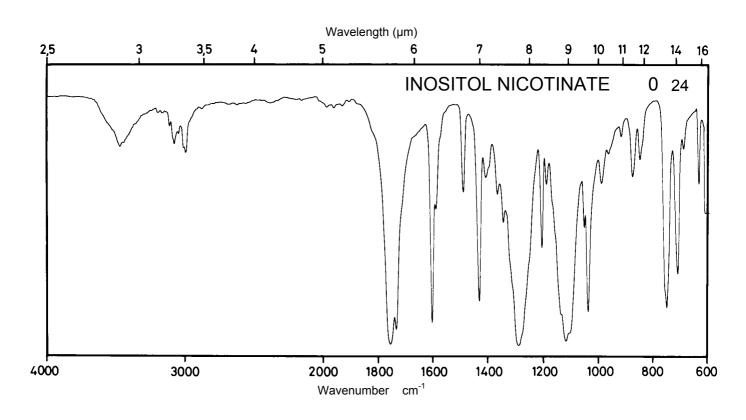



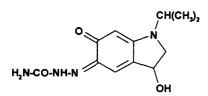



264.3

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	357 nm 221 nm	356 nm 221 nm	378 nm	448 nm
E ^{1%} 1cm	1080 672	1014 515	1200	1335
3	28600 17800	26800 13600	31800	35300

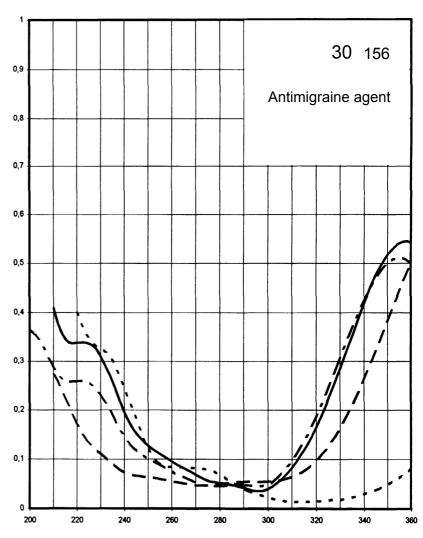

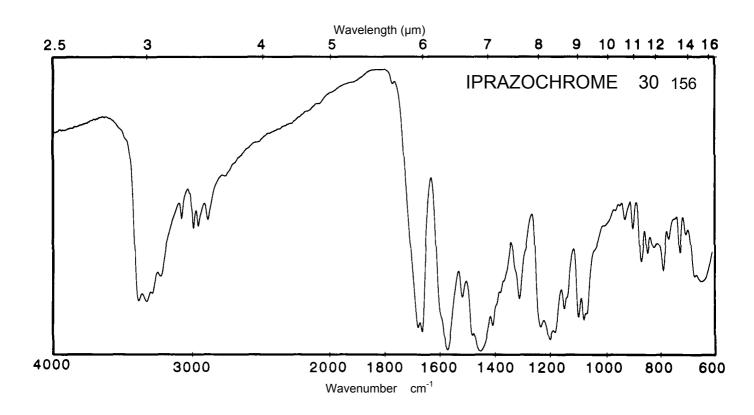




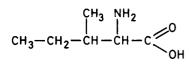



131.2


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

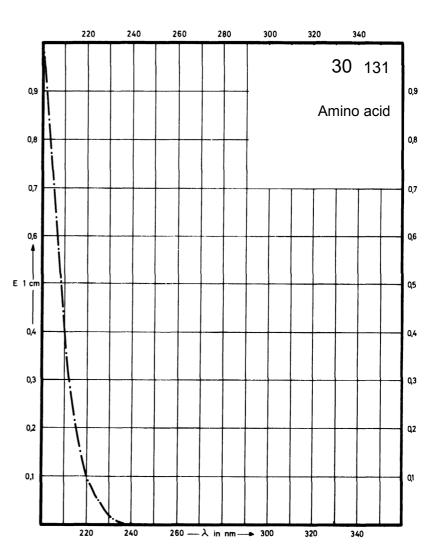


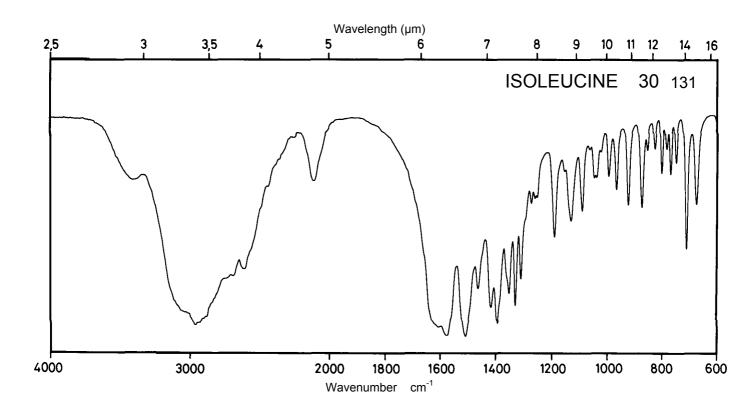

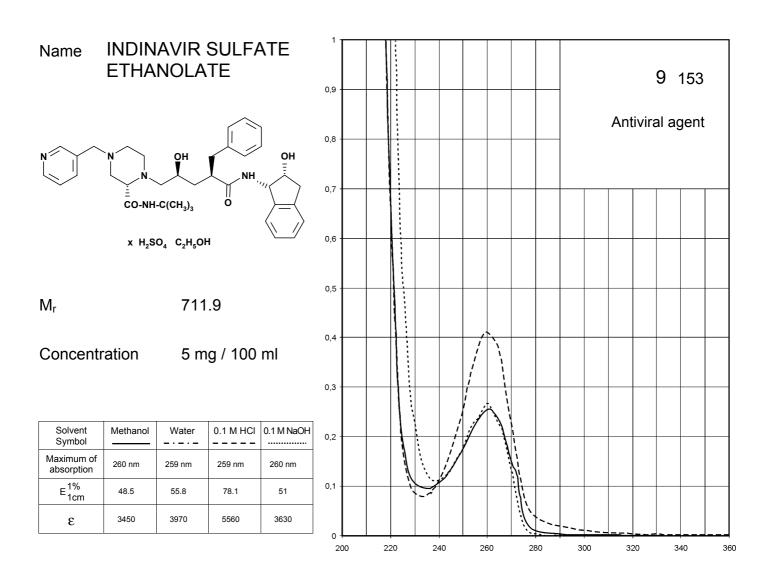

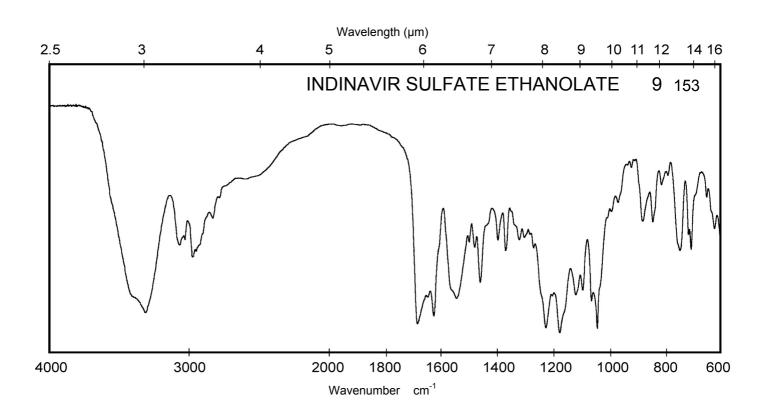



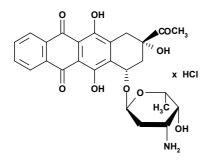



 $\,M_r\,$

534.0


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	481 nm 287 nm 251 nm	484 nm 289 nm 257 nm	484 nm 289 nm 257 nm	Decom- position observed
E ^{1%} 1cm	207 179 816	194 180 743	194 180 738	
ε	11100 9540 43600	10400 9620 39700	10400 9620 39400	

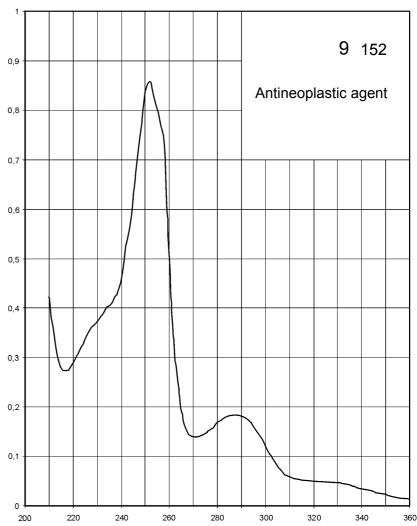

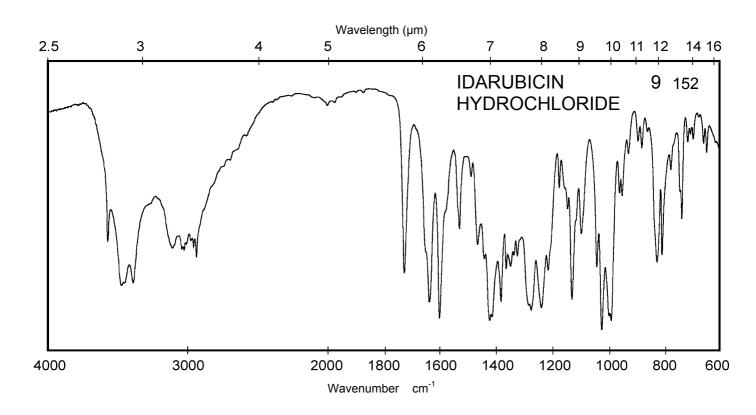




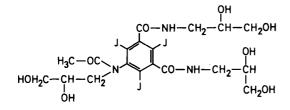


821.1


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm		244 nm	243 nm
E ^{1%} 1cm	359		355	357
3	29500		29200	29300

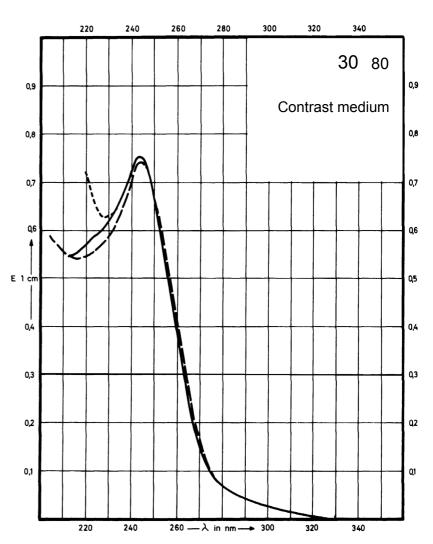

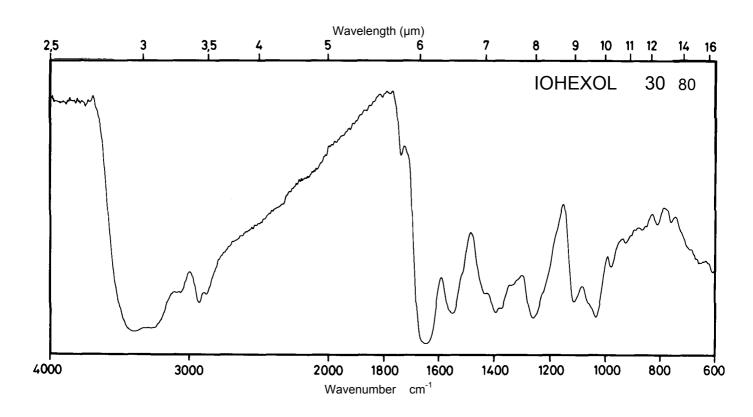




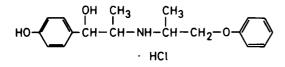


 M_r 337.8

Concentration 2 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 270 nm 221 nm		276 nm 270 nm	292 nm 277 nm 244 nm
E ^{1%} 1cm	83 78 436		71 71	70 71 396
3	2800 2620 14720		2400 2400	2350 2410 13380

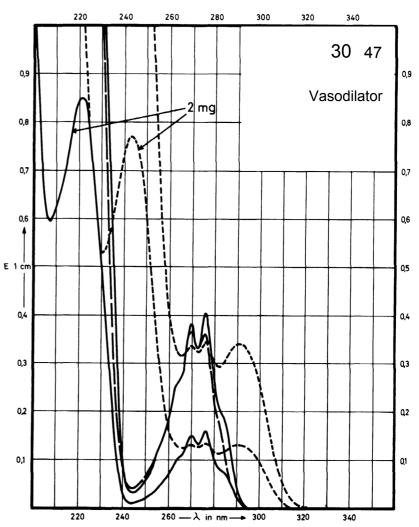

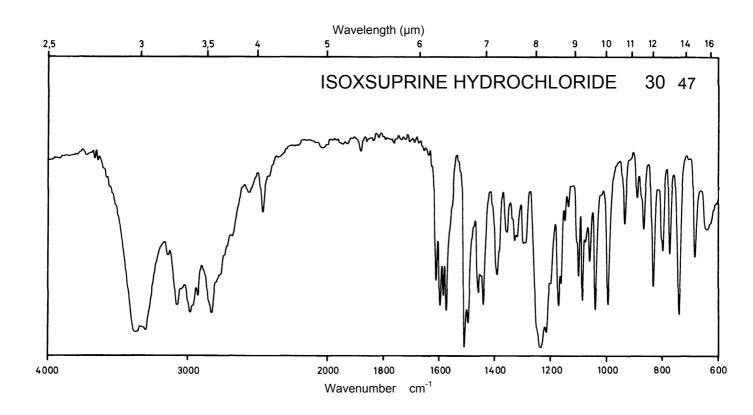

Name ISONICOTINIC ACID


Соон


 M_{r}

123.1


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		270 nm	266 nm
E ^{1%} 1cm	215		362	199
3	2650		4460	2450

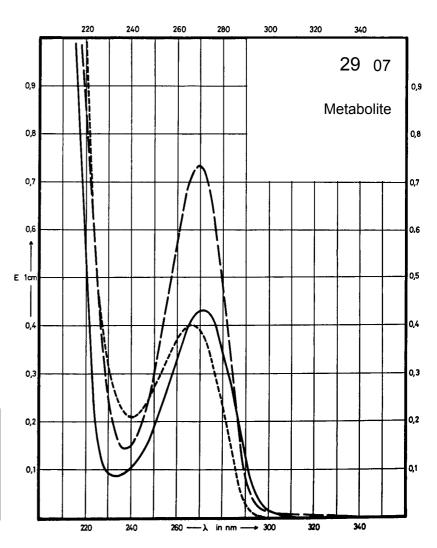

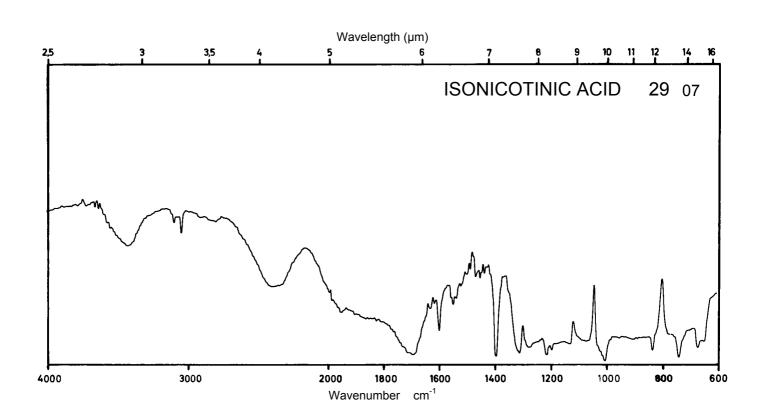




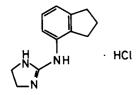


237.7


Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm 266 nm		272 nm 265 nm	
E ^{1%} 1cm	43 47		47 51	
ε	1020 1110		1110 1200	

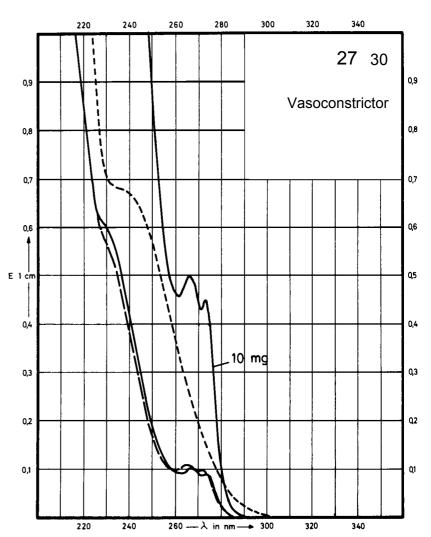

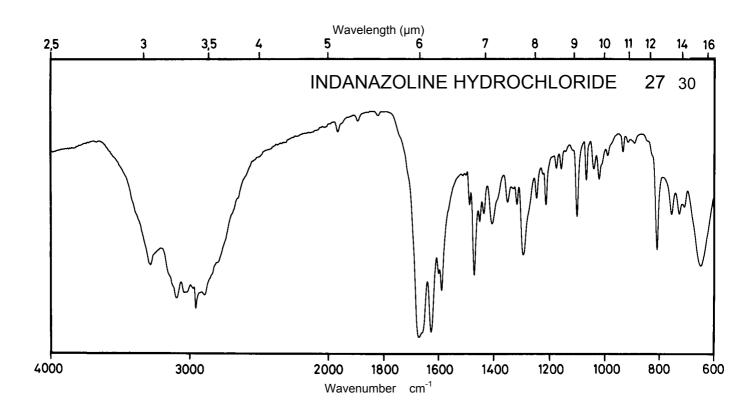




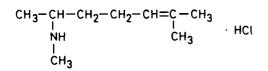


 M_{r} 177.7


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

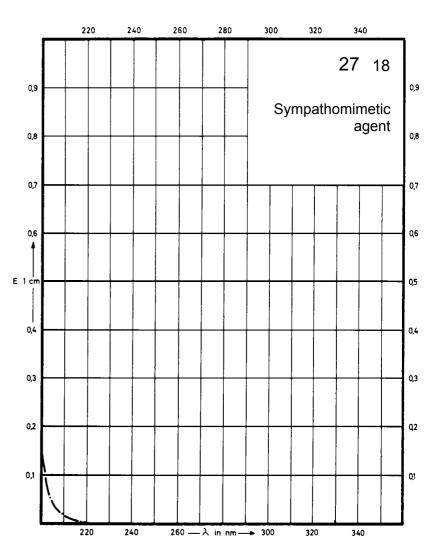

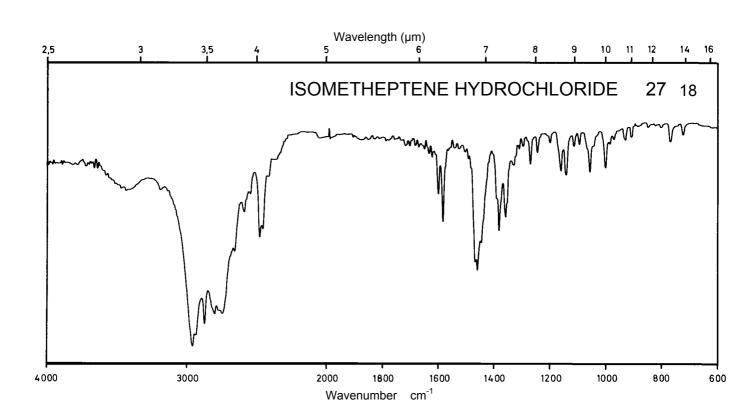

Name ISOPRENALINE HYDROCHLORIDE

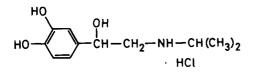



 M_r

247.7


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm 225 nm		279 nm 220 nm	
E 1% 1cm	129 239		112 253	
ε	3200 5900		2770 6270	

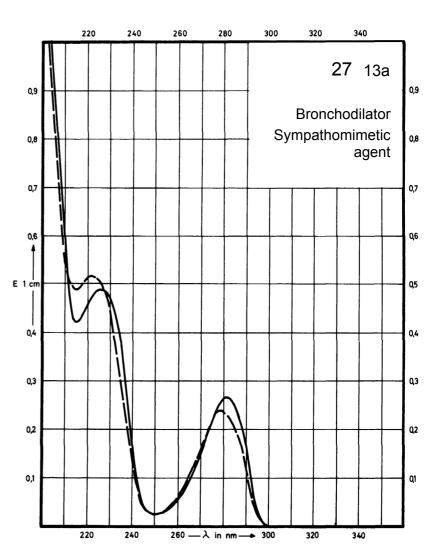

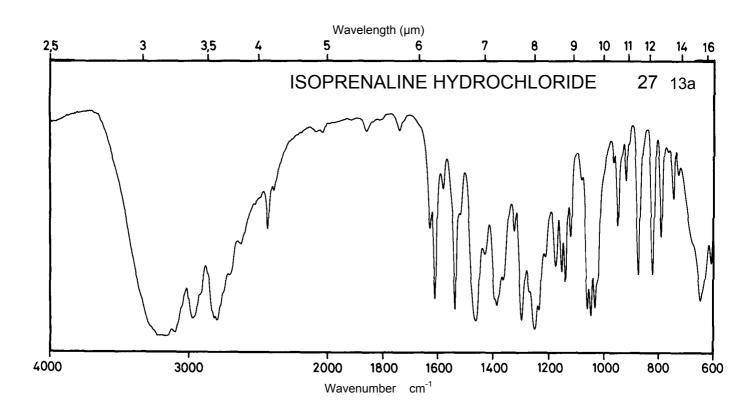




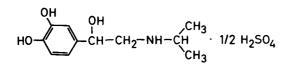


260.3


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	f 282 nm		279 nm	Decom- position observed
E ^{1%} _{1cm}	121		106	
ε	3150		2760	

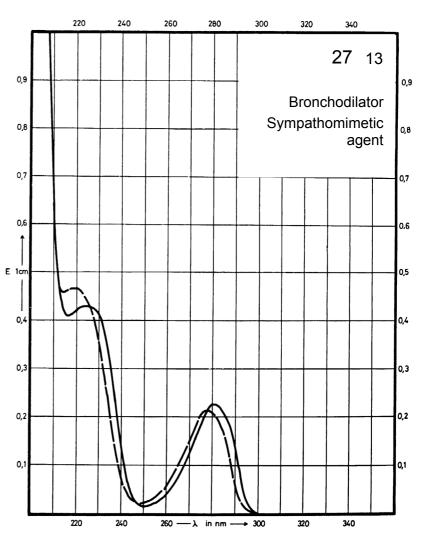

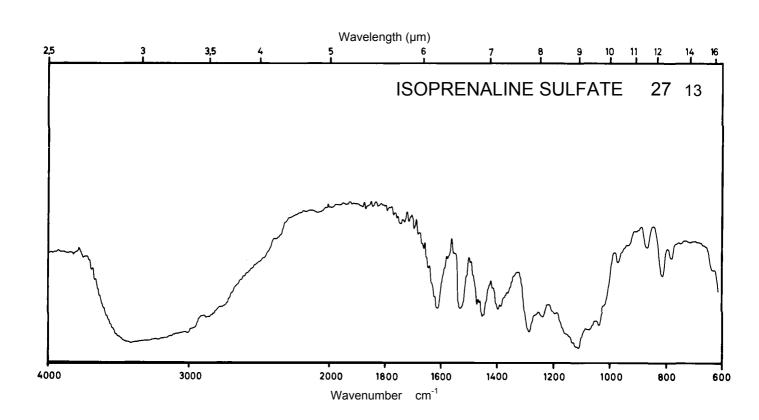




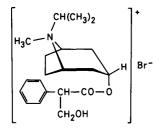


412.4


Concentration 10 mg / 100 ml 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm	263 nm 257 nm 252 nm		
E ^{1%} 1cm	3.30 4.36 3.56	3.42 4.52 3.74		
3	136 180 147	140 187 154		

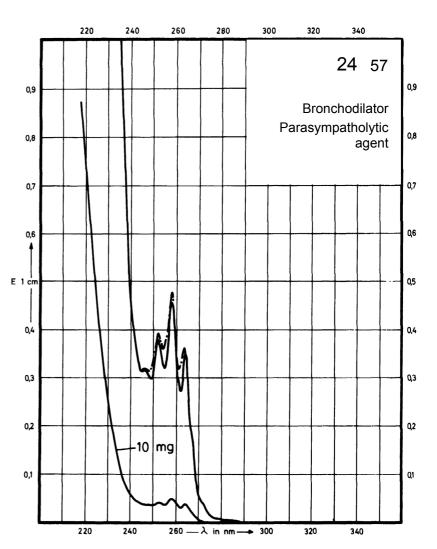

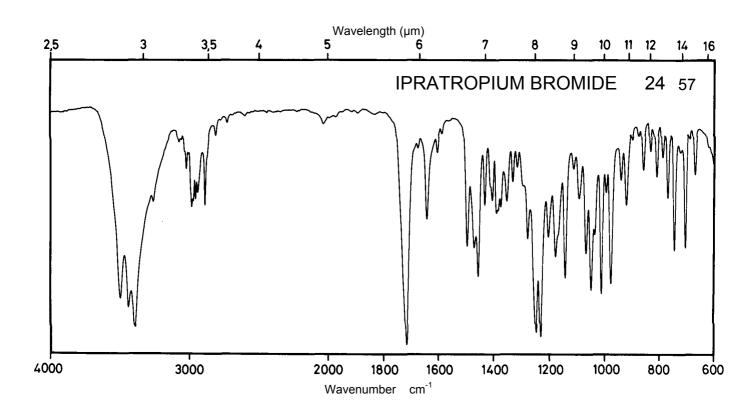




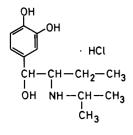


275.8


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm		280 nm	
E ^{1%} 1cm	117		104	
ε	3220		2860	

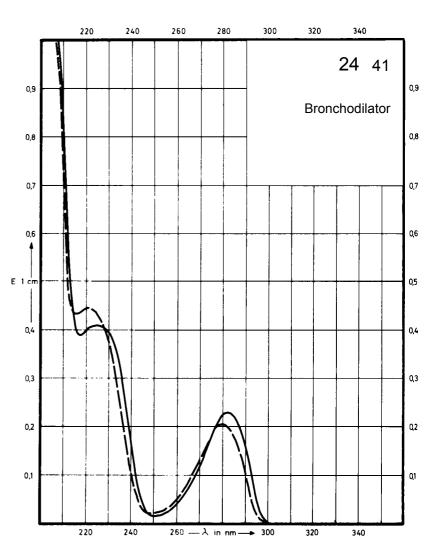

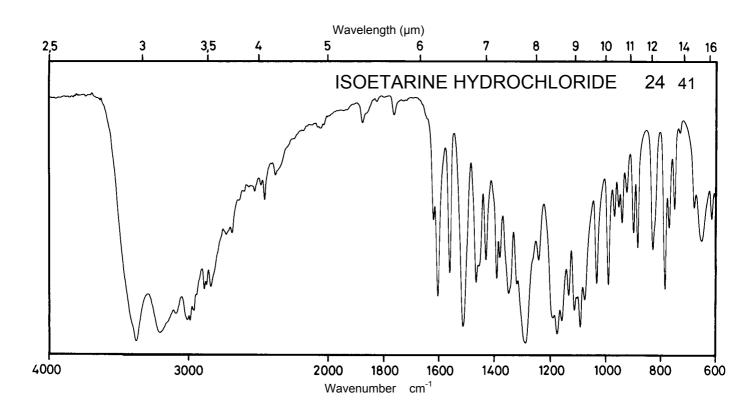




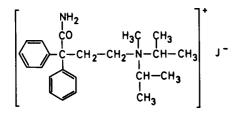



 M_r 480.4


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm 260 nm			
E ^{1%} 1cm	6.8 9.1			
3	326 436			

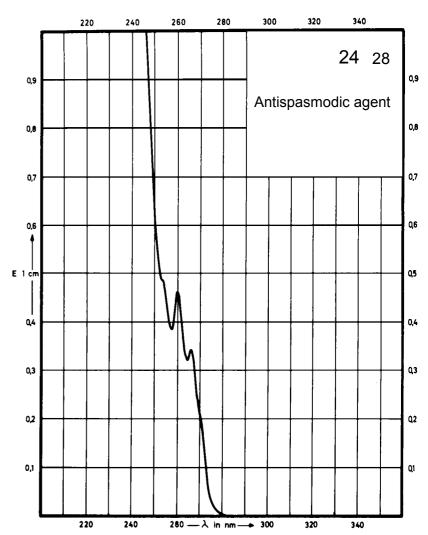

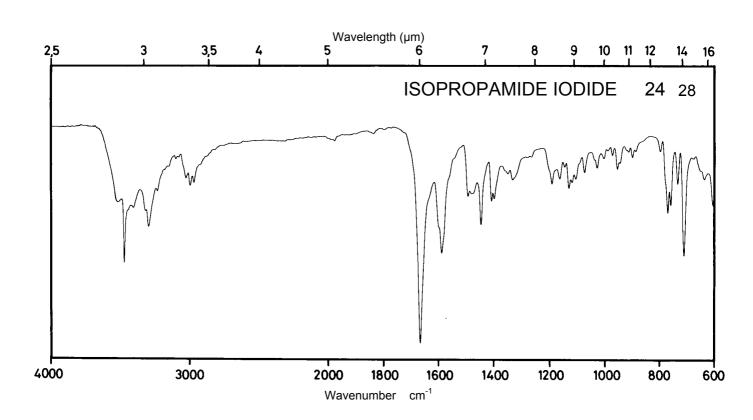

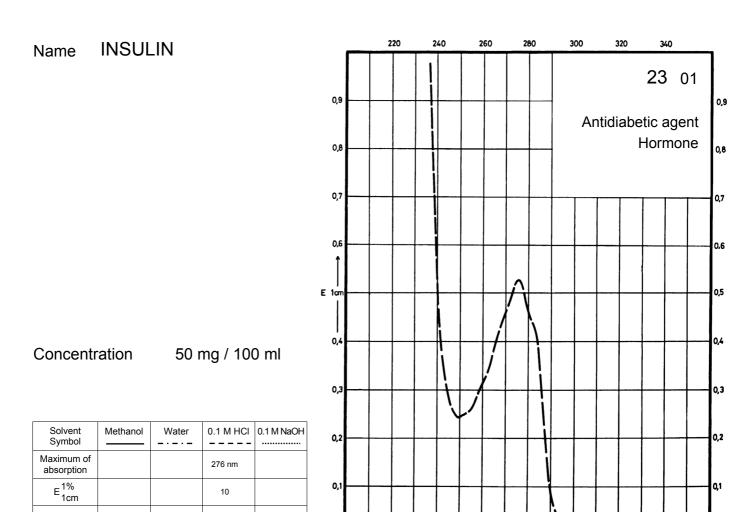

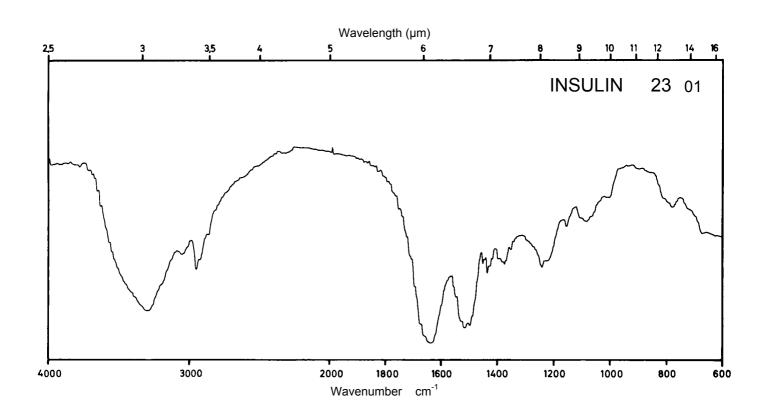

−A in nm

Name IPRONIAZID PHOSPHATE


 M_{r}


277.2

Concentration


2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm		268 nm	Decom- position observed
E ^{1%} 1cm	169		184	
3	4680		5091	

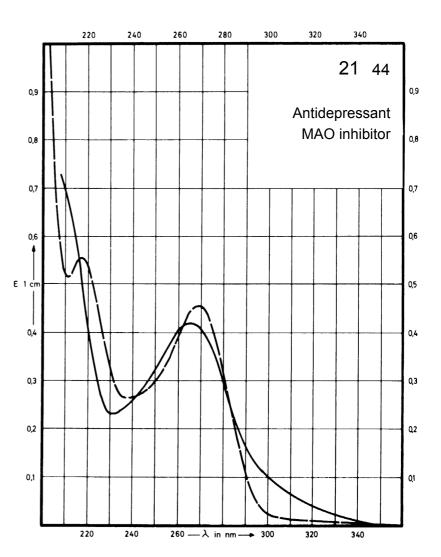

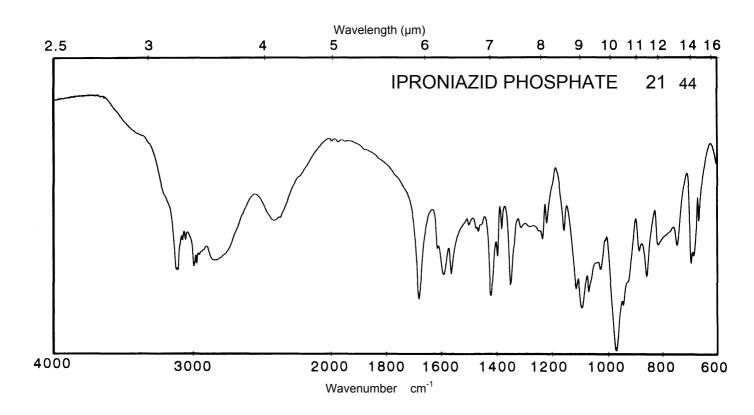

Name IMIPRAMINE HYDROCHLORIDE

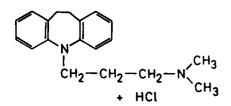



 M_{r}

316.9


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	251 nm		249 nm	252 nm
E ^{1%} 1cm	274		255	224
ε	8680		8080	7100

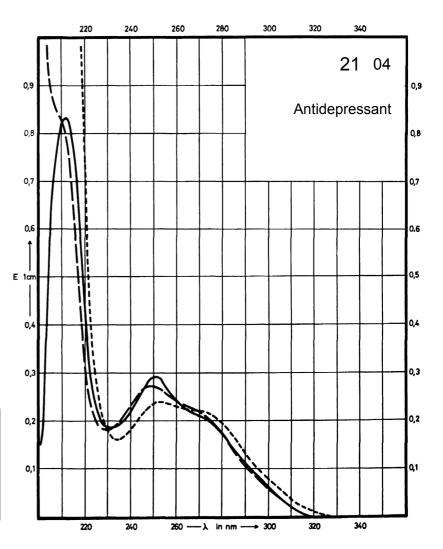

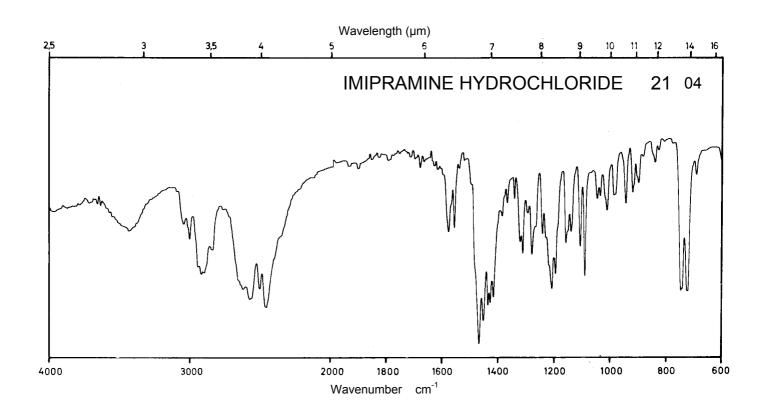




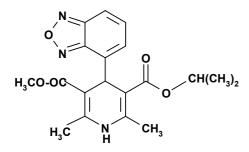


371.4


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	326 nm 228 nm	332 nm 239 nm	331 nm 239 nm	332 nm
E ^{1%} 1cm	286 470	266 442	262 435	262
3	10600 17400	9860 16400	9720 16200	9750

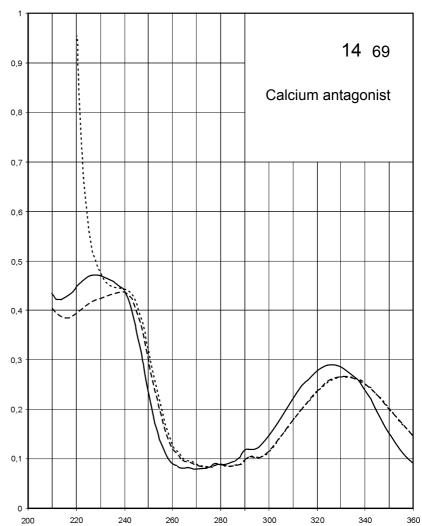

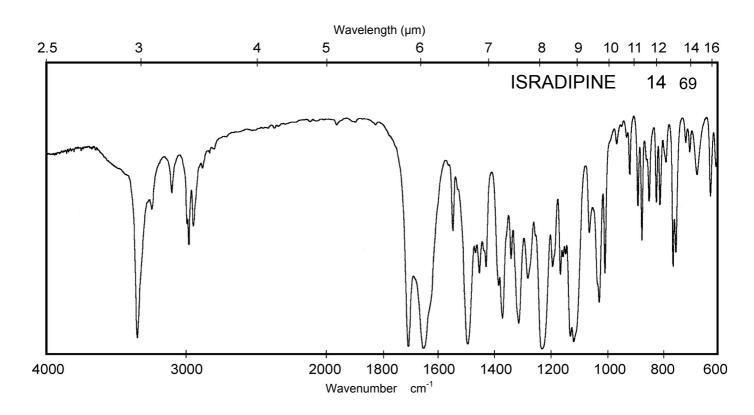



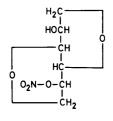



 M_r

191.1


Concentration 6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

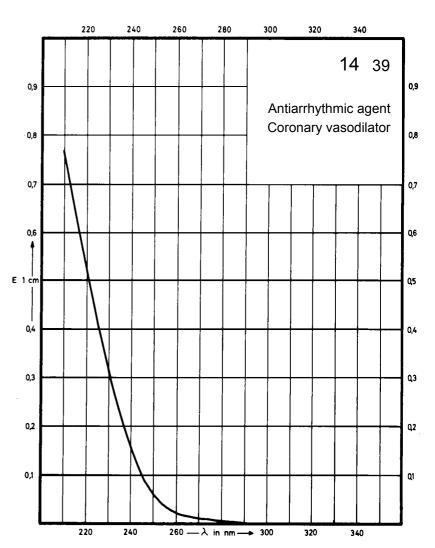

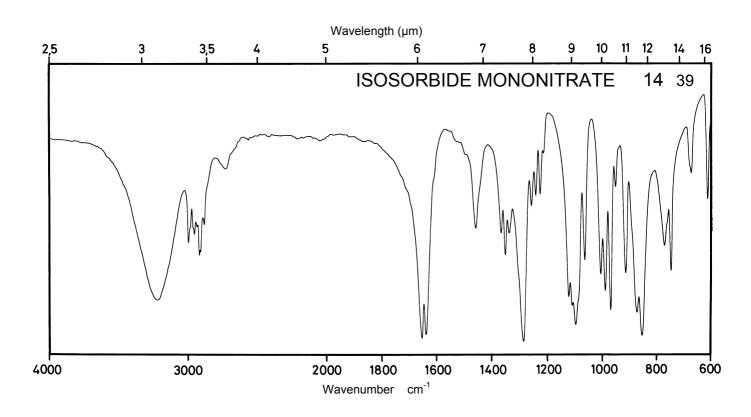



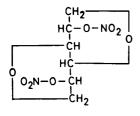



 M_r

236.1


Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

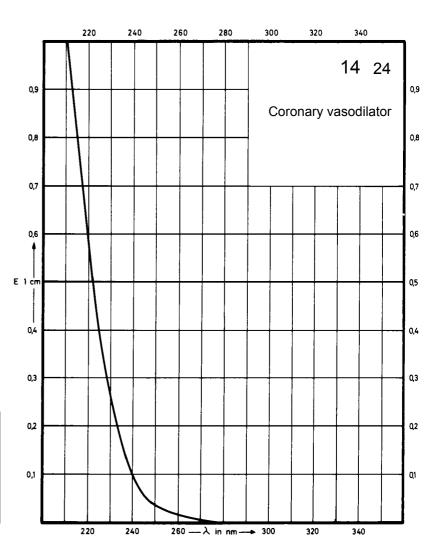

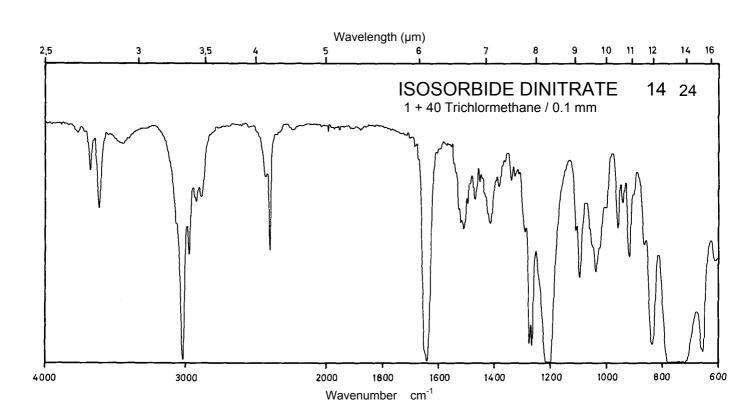




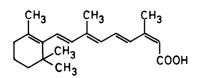



300.4

Concentration

0.31 mg / 100 ml

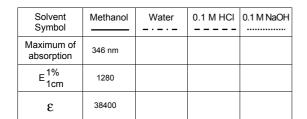

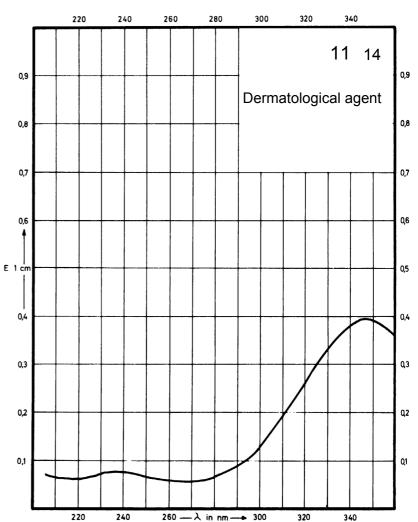

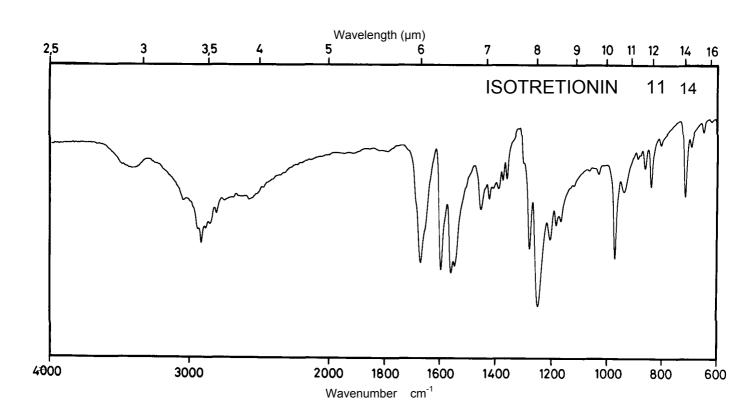


Name ISOCONAZOLE NITRATE


 M_{r}

479.1

Concentration 52 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm 273 nm			
E ^{1%} 1cm	13.6 15.5			
3	650 740			

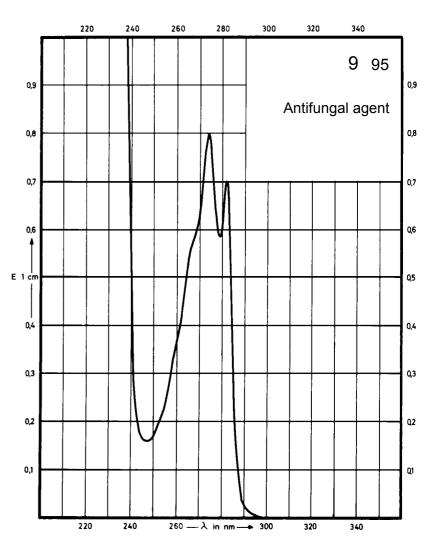

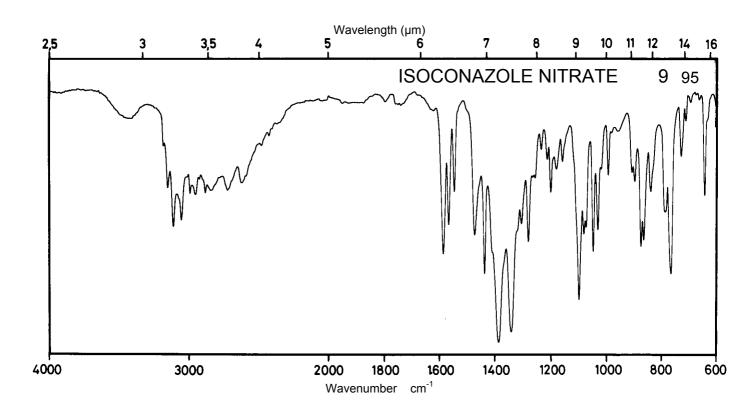


Name ISOCONAZOLE


 M_{r}

416.1


Concentration 50 mg / 100 ml

	,			
Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 273 nm			
E ^{1%} 1cm	15.3 17.5			
ε	640 730			

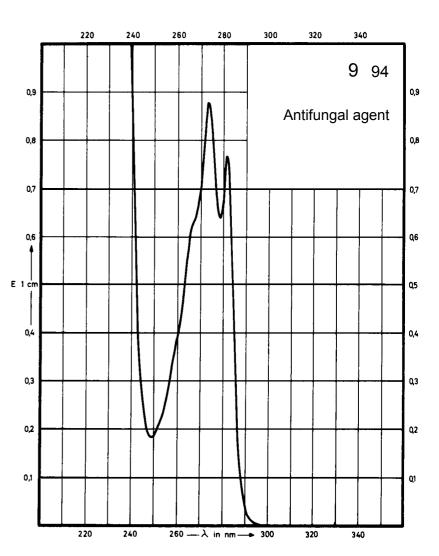

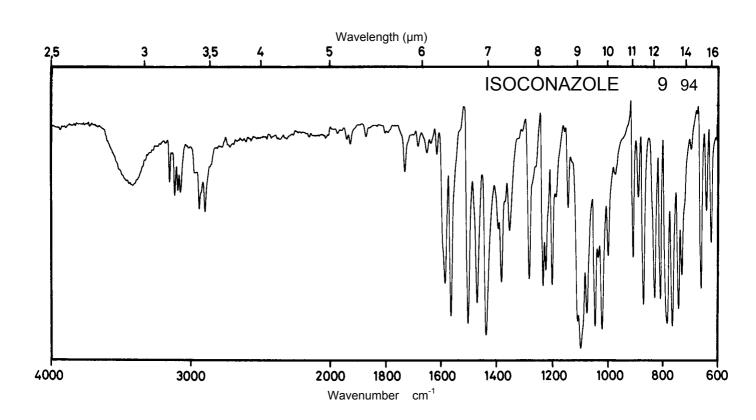




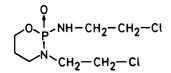


261.1

Concentration 60 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
3				

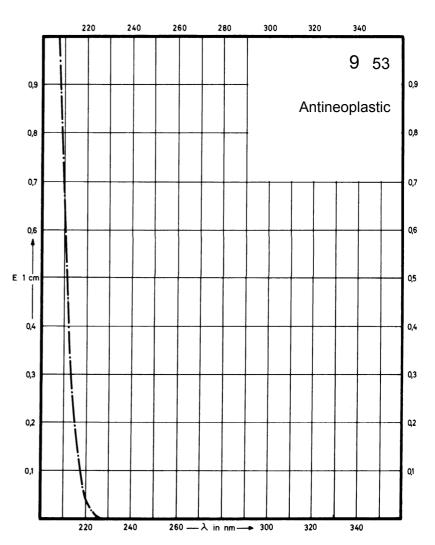

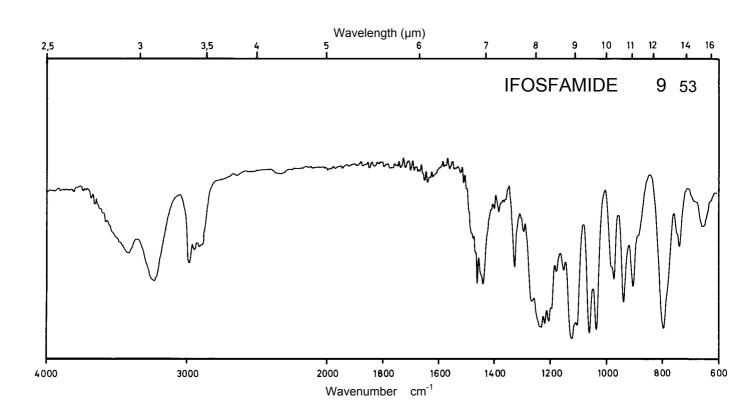


Name IDOXURIDINE


0 = N - ÇH - CH2 CH2 CH0H CH2OH

 M_{r}


354.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm	288 nm	288 nm	280 nm
E ^{1%} 1cm	227	217	217	164
ε	8040	7680	7680	5810

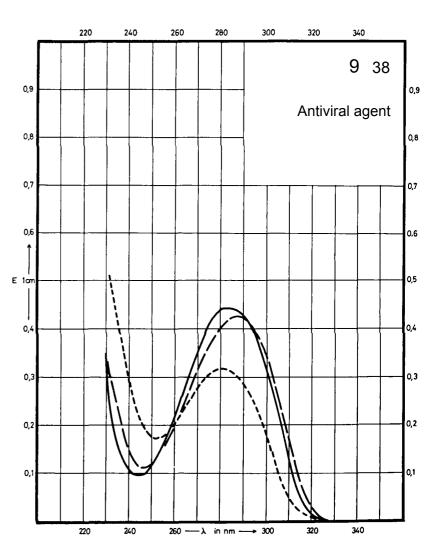

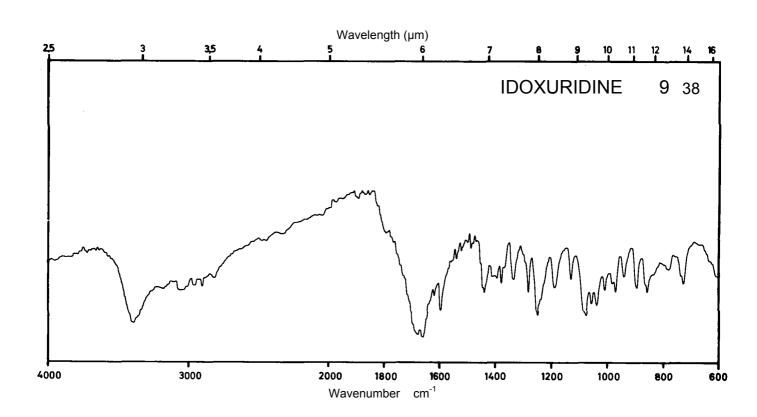





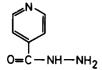


137.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Cyllibol				
Maximum of absorption	262 nm		266 nm	296 nm
E ^{1%} 1cm	334		380	368
3	4580		5210	5050

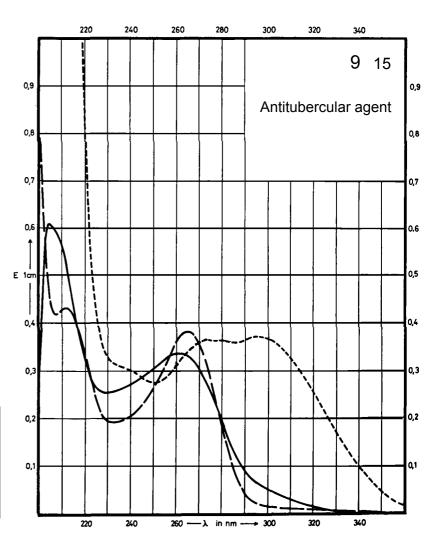

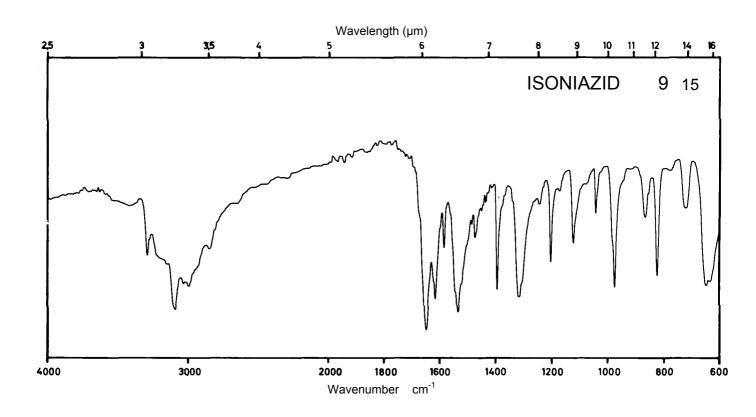




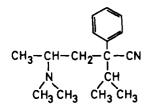


 M_{r}


244.4

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	265 nm		264 nm	
	259 nm		258 nm	
absorption	253 nm		252 nm	
10/_	6.8		6.2	
E ^{1%} 1cm	8.6		7.8	
1cm	7.6		6.0	
	166		152	
3	209		191	
	195		1/17	

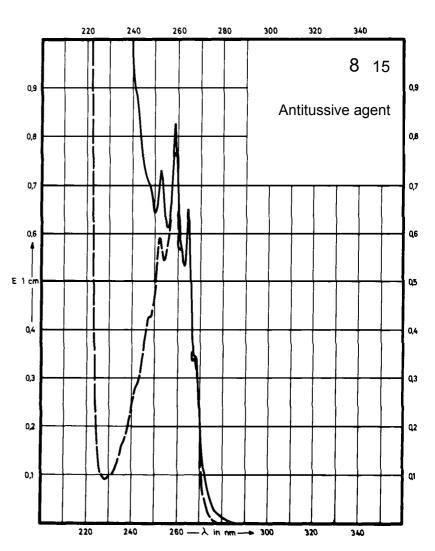

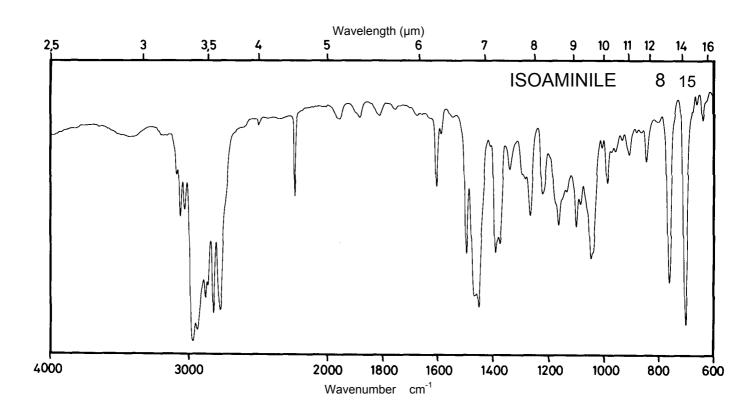




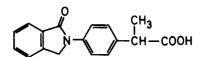


 M_r 281.3


Concentration 1 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 281 nm 227 nm 276 nm 227 nm 279 nm absorption E^{1%} 1cm 478 486 515 547 456 14500 15400 13400 13700 12800 3

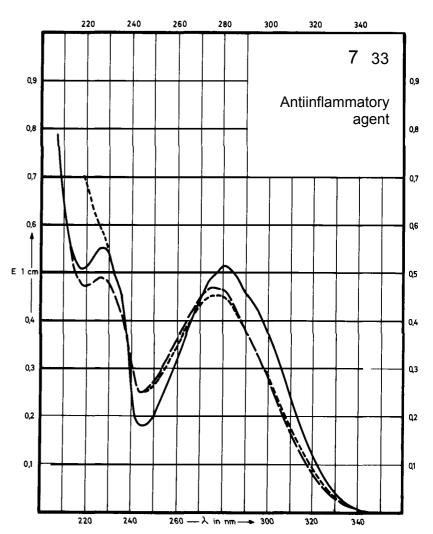

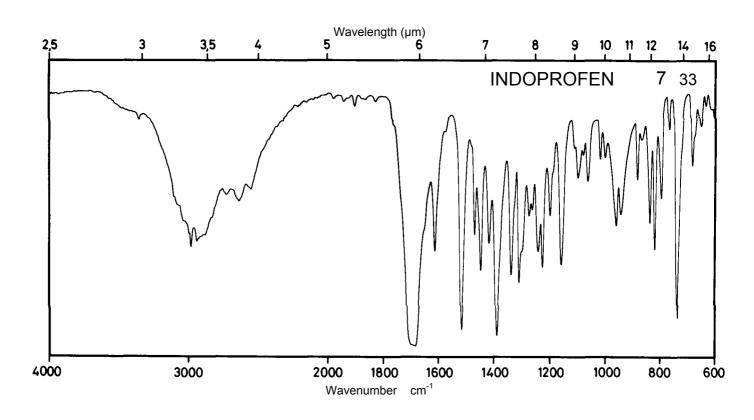





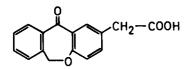


 M_r 268.2

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	336 nm		337 nm	342 nm
	270 nm		275 nm	276 nm
absorption	224 nm		222 nm	2/011111
1%	117		123	120
E ^{1%}	552		547	560
TCITI	806		760	300
	3130		3300	3210
3	14800		14670	15030
	21620		20380	15050

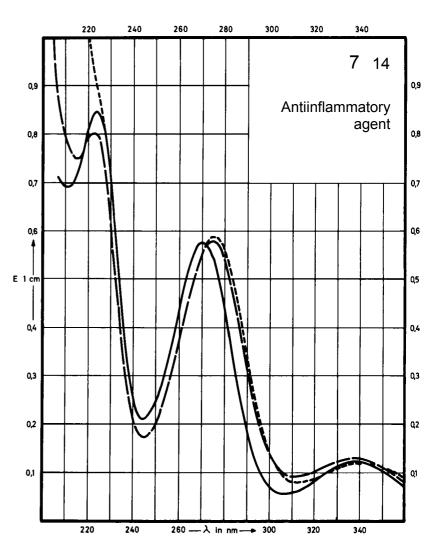

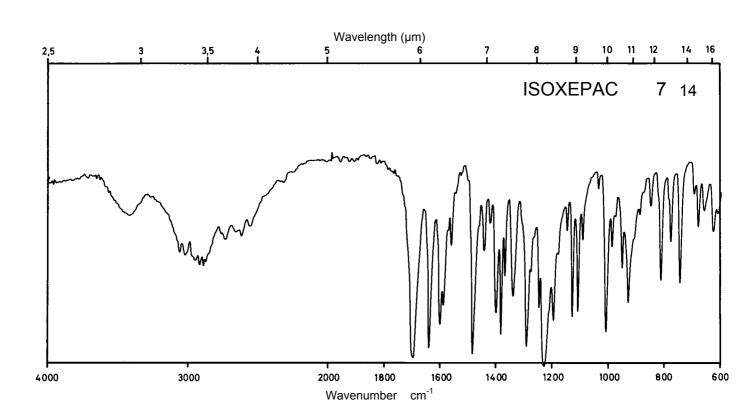





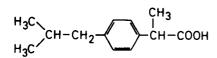

206.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 264 nm 258 nm			272 nm 264 nm 258 nm
E ^{1%} 1cm	11.2 14.5 11.3			15.4 18.4 15.0
ε	230 300 233			320 380 310

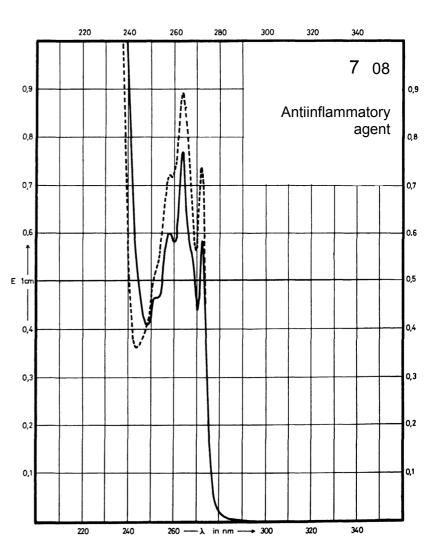

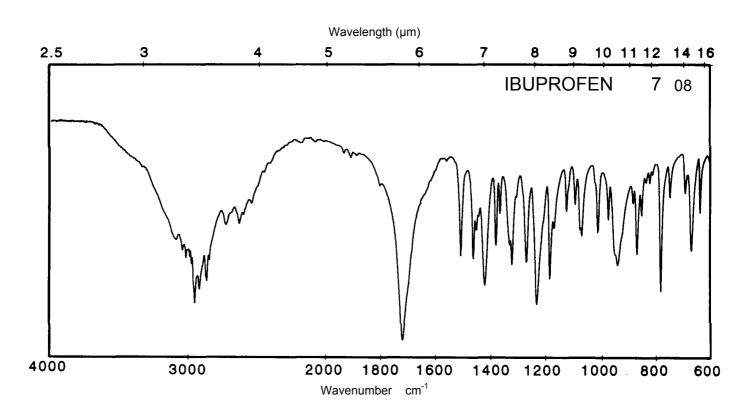



Name INDOMETACIN

 M_{r}


357.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm		316 nm 260 nm	280 nm
E ^{1%} 1cm	172		177 487	201
3	6150		6330 17420	7190

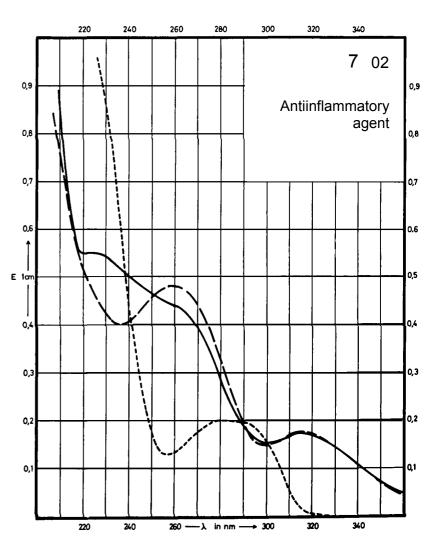

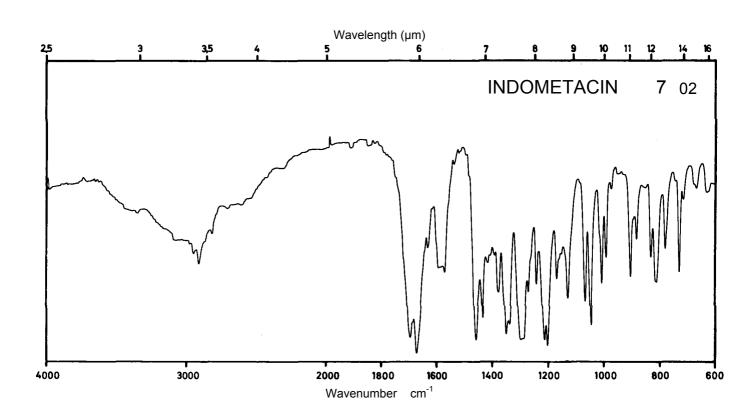



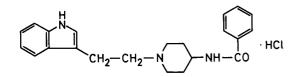


383.9


2.5 mg / 100 ml


 M_r


Concentration

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 287 nm 277 nm 271 nm 125 157 289 nm 279 nm 273 nm Maximum of 288 nm 278 nm absorption 142 170 E^{1%} 1cm 124 149 168 5460 4760 5720 3 6520 6020 6020 6450

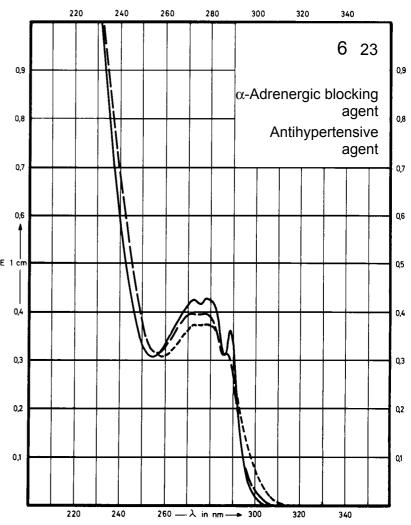


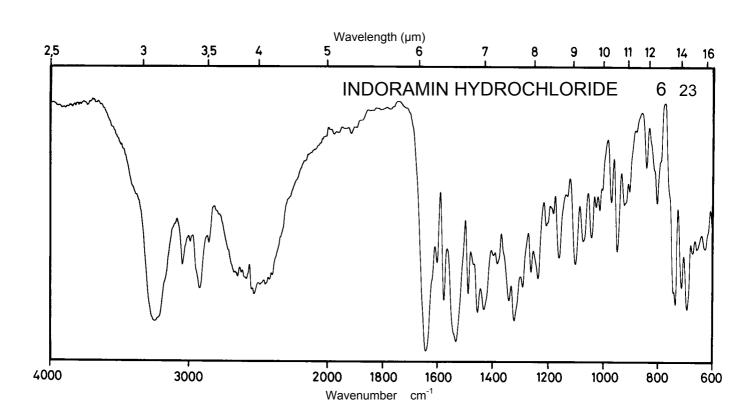

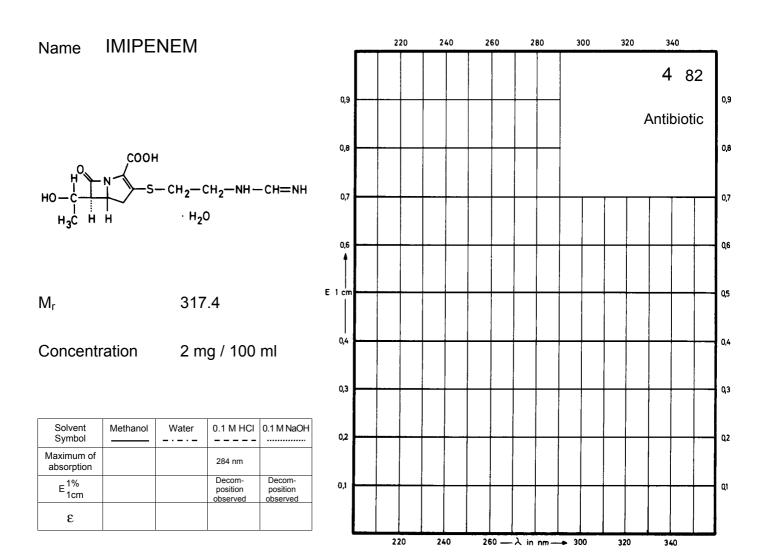

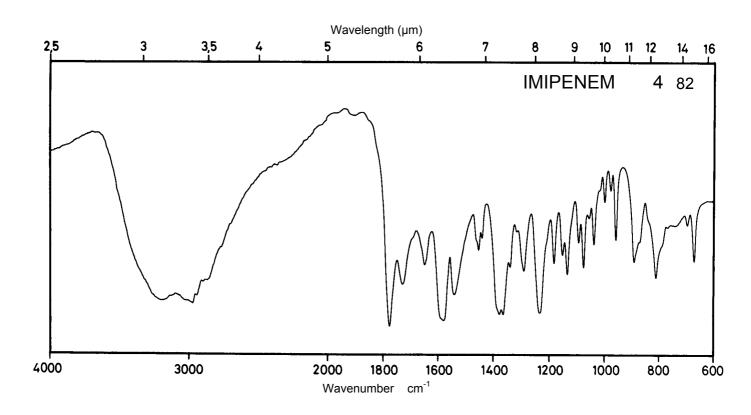


Name ISOTHIPENDYL HYDROCHLORIDE

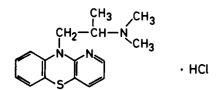


 M_{r}


321.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 248 nm		330 nm 244 nm	320 nm 250 nm
E ^{1%} 1cm	136 873		112 756	145 776
ε	4380 28100		3610 24340	4670 24980

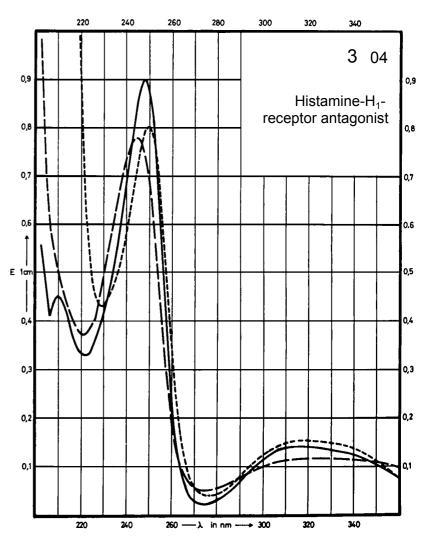

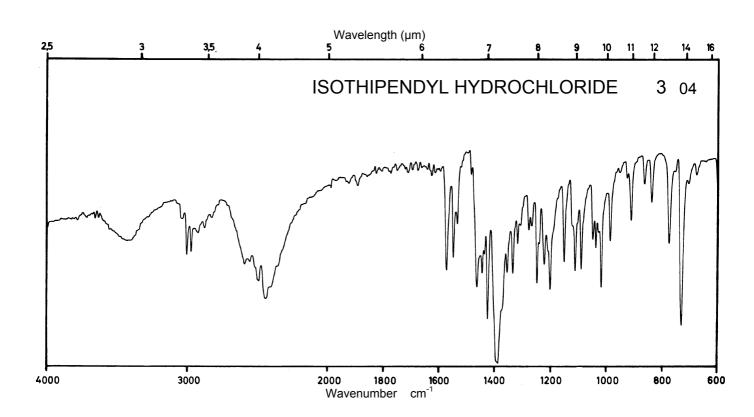





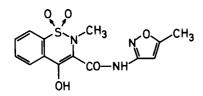



335.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				347 nm 248 nm
E ^{1%} 1cm				393 440
3				13200 14800

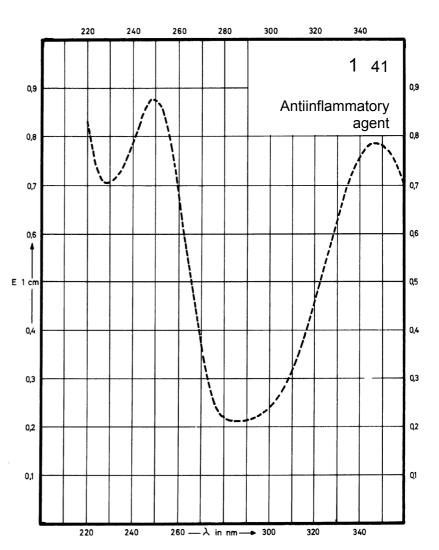


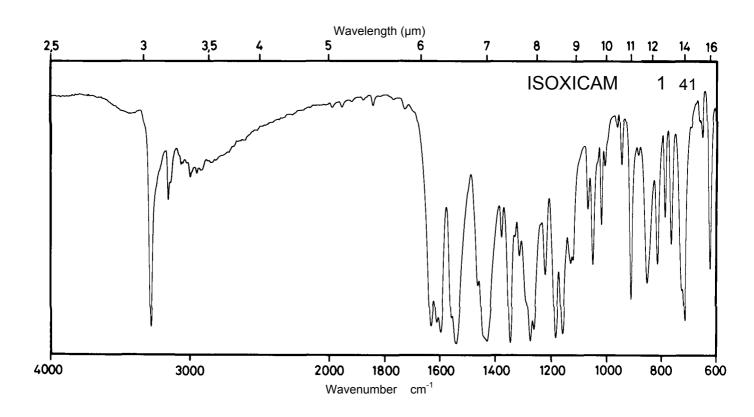

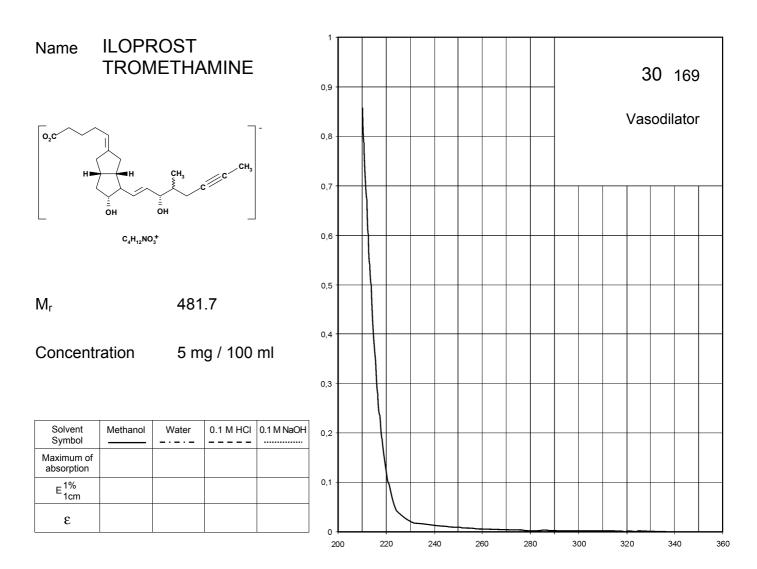

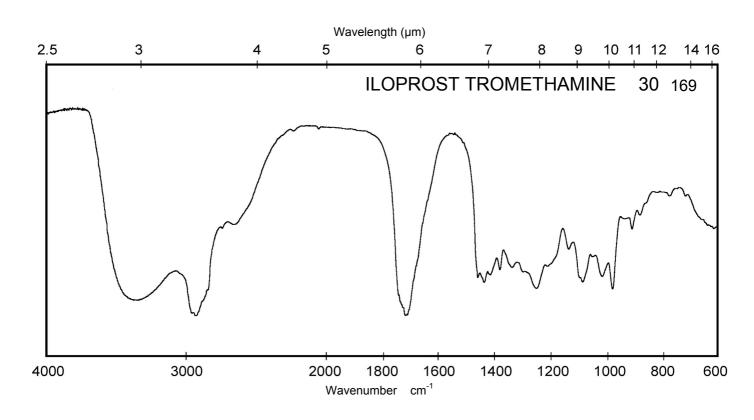




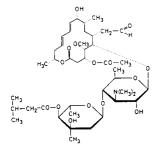


828.0


Concentration 2.6 mg / 100 ml 52 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			284 nm 231 nm	
E ^{1%} 1cm			4.3 312	
ε			350 25850	

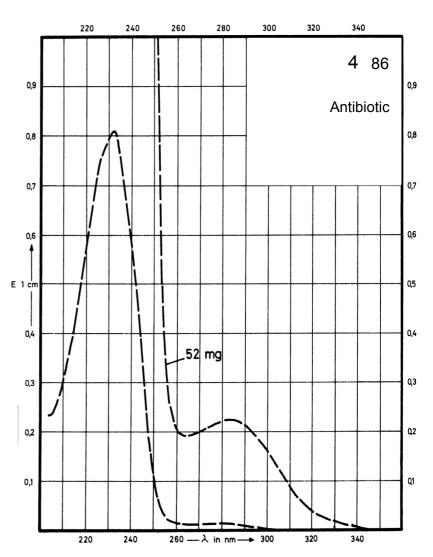

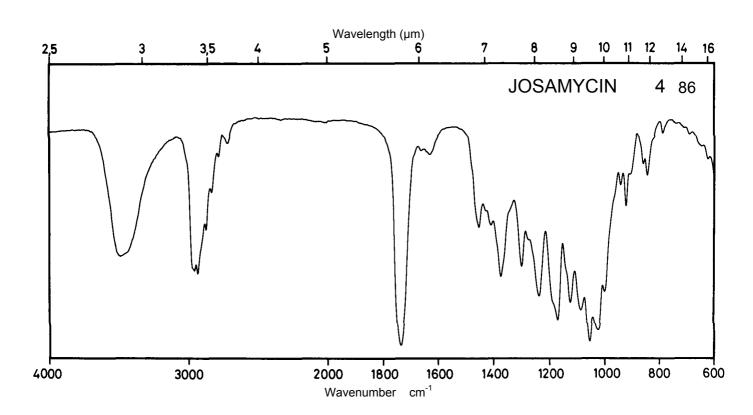




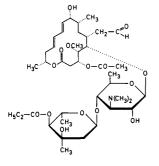


 $\,M_{r}\,$

0.008

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		231 nm		
E ^{1%} _{1cm}		311		
3		24900		

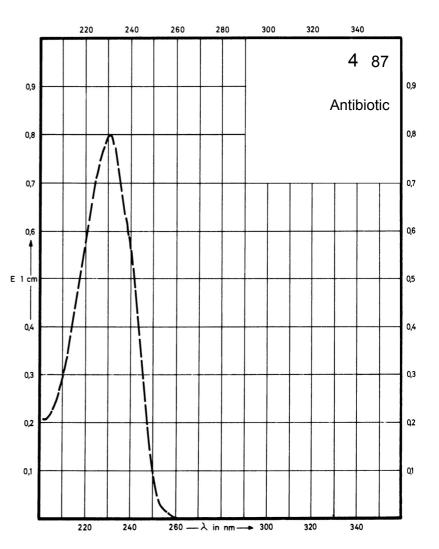

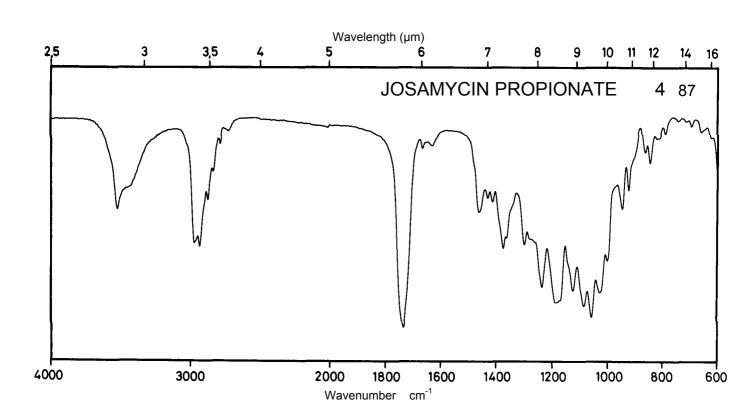



Name KETOBEMIDONE HYDROCHLORIDE

 M_{r}


283.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm		280 nm	300 nm
E ^{1%} 1cm	87		76	110
ε	2460		2160	3110

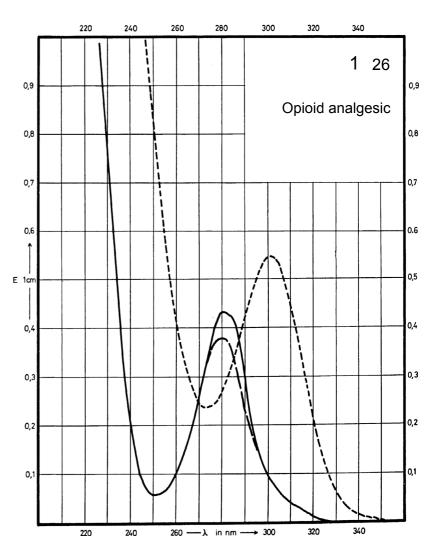

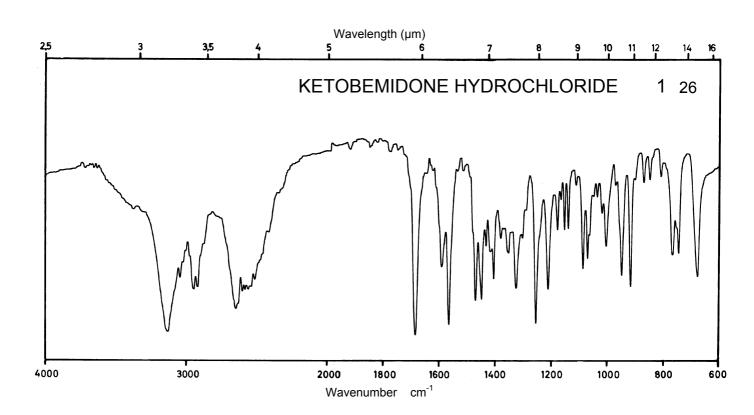





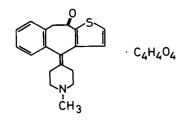


 M_r 425.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	297 nm	301 nm	302 nm	305 nm
E ^{1%} 1cm	325	334	338	298
ε	13800	14230	14400	12680

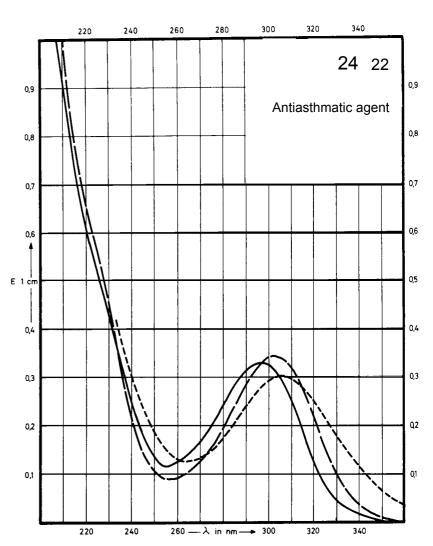

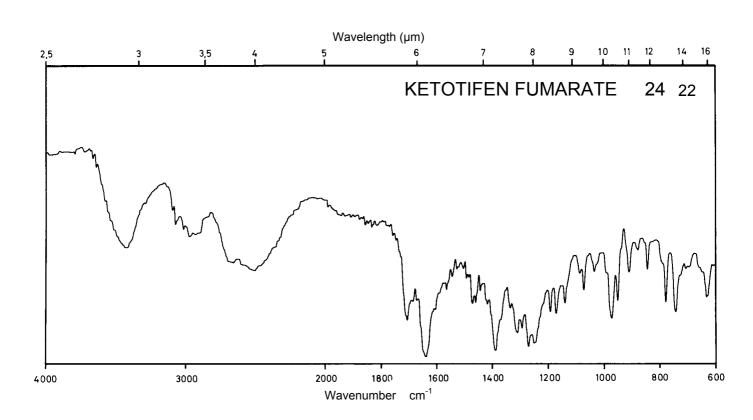




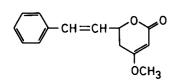


230.3


Concentration 0.5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	292 nm 283 nm 245 nm		247 nm	Decom- position observed
E ^{1%} 1cm	46 67 1190		1300	
ε	1070 1550 27400		29940	

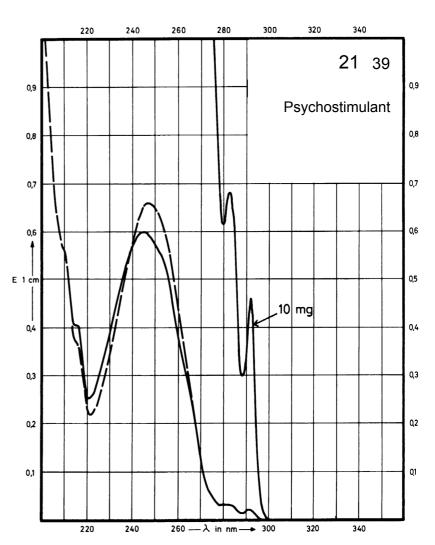

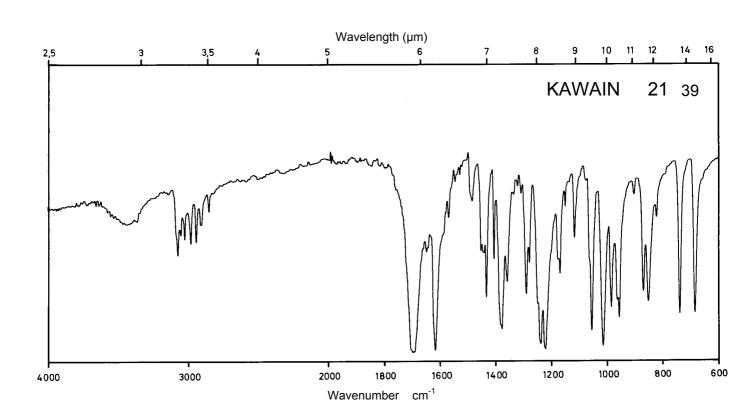




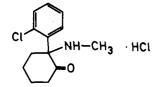


Concentration 20 mg / 100 ml


274.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm	277 nm	277 nm	297 nm
	270 nm	270 nm	270 nm	268 nm
E ^{1%} 1cm	19.7	20.2	20.2	4.6
	22.3	23.1	23.1	9.9
3	540	550	550	126
	610	630	630	270

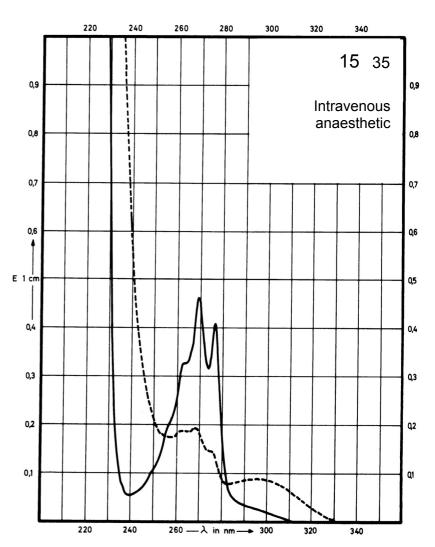

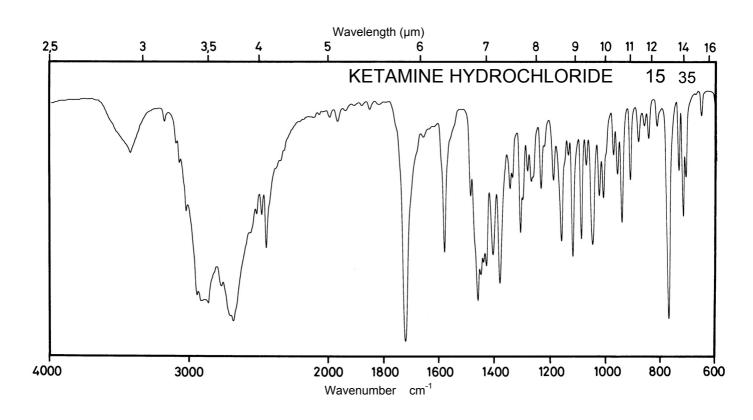




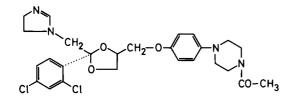


 M_r 531.4


Concentration 1.25 mg / 100 ml 2.5 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm 244 nm		270 nm 223 nm	287 nm
E ^{1%} 1cm	32 280		27 530	29
3	1700 14900		1420 28000	1500

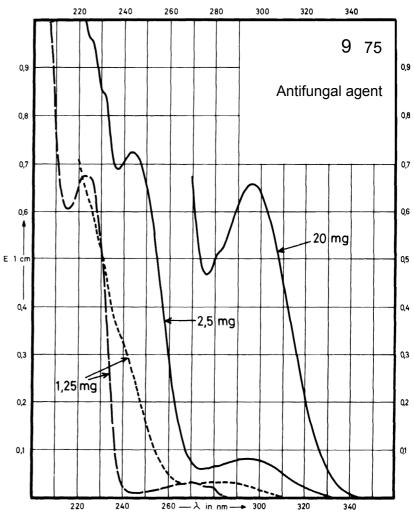

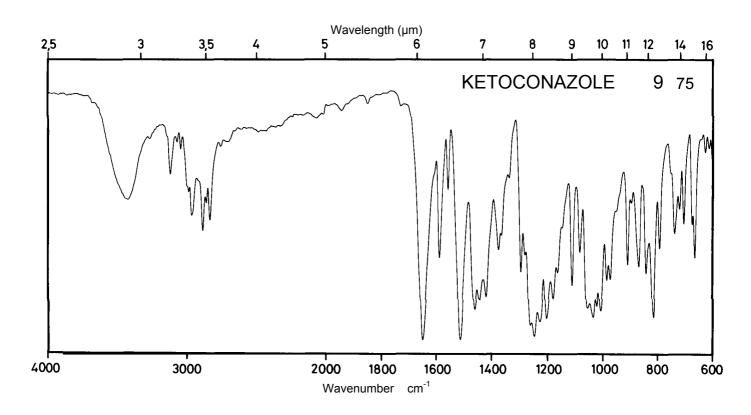





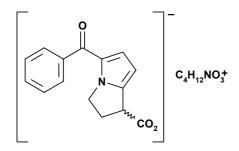


376.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	319 nm	323 nm	318 nm	323 nm
	244 nm	250 nm	250 nm	249 nm
E ^{1%}	531	544	510	544
1cm	203	184	184	184
ε	20000	20500	19200	20500
	7660	6940	6940	6940

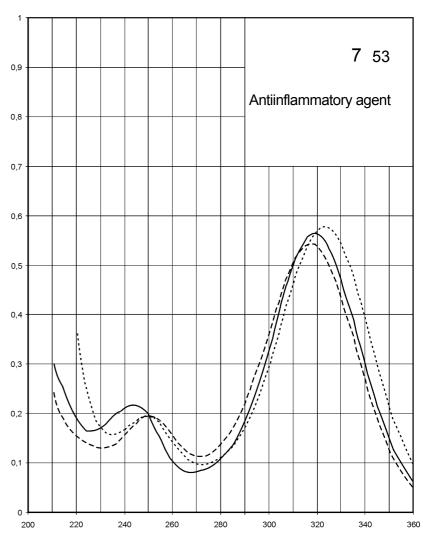

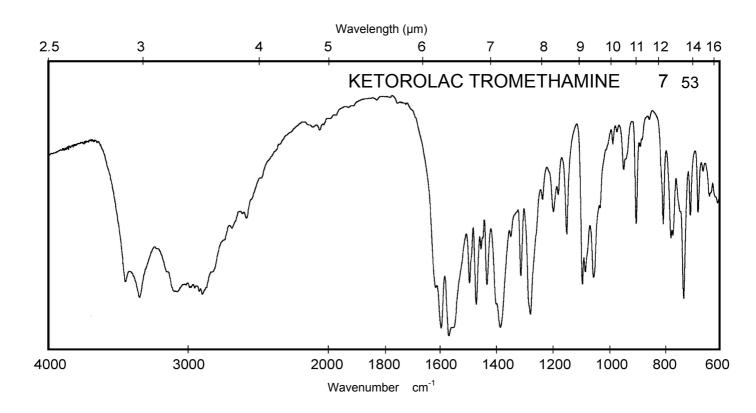





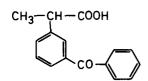

254.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	255 nm		260 nm	261 nm
E ^{1%} 1cm	672		666	649
ε	17080		16930	16500

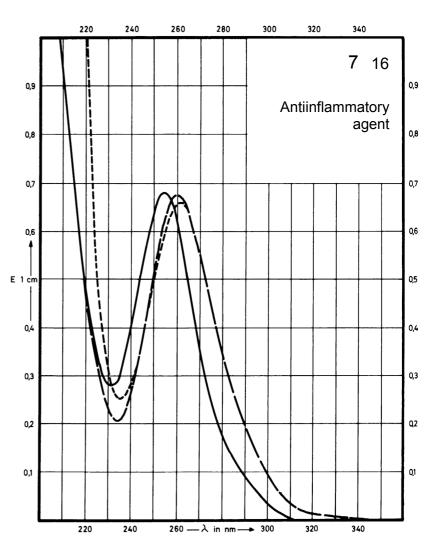

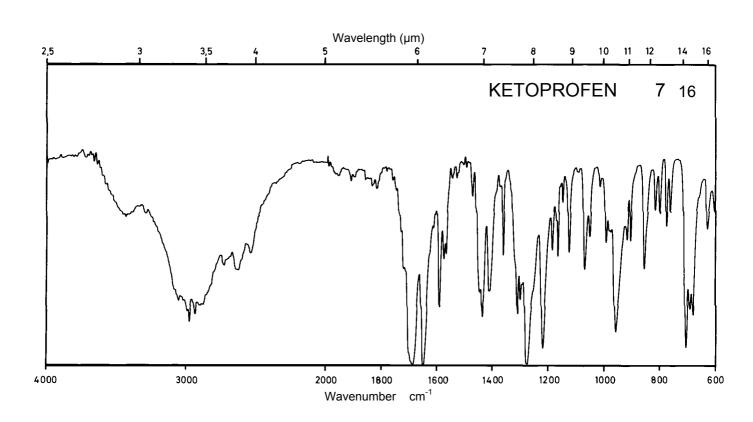



Name KEBUZONE

 M_r

322.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	244 nm		237 nm	262 nm
E ^{1%} 1cm	448		404	617
3	14440		13020	19890

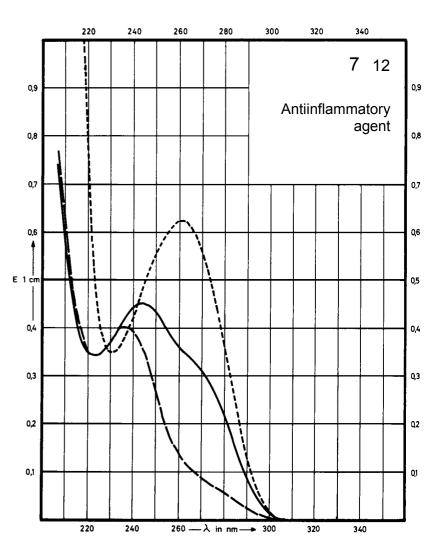

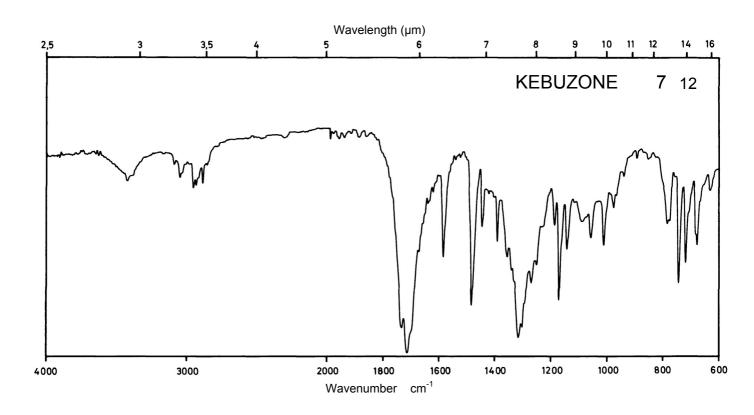

Name KHELLIN

H₃C-0 0 CH₃


 M_{r}


260.3


Concentration

0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	330 nm 247 nm		335 nm 249 nm	335 nm 249 nm
E ^{1%} 1cm	180 1470		175 1505	175 1485
ε	4690 38260		4560 39180	4560 38650

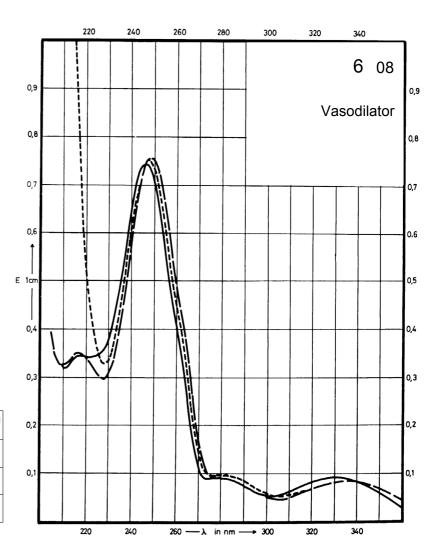


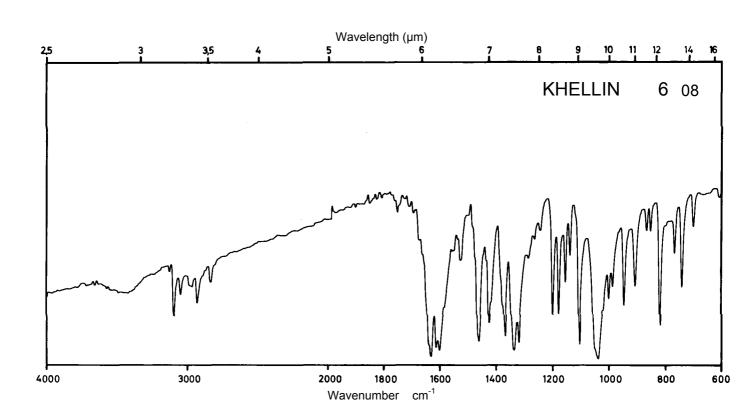


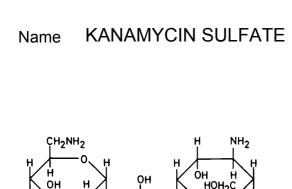


582.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
3				

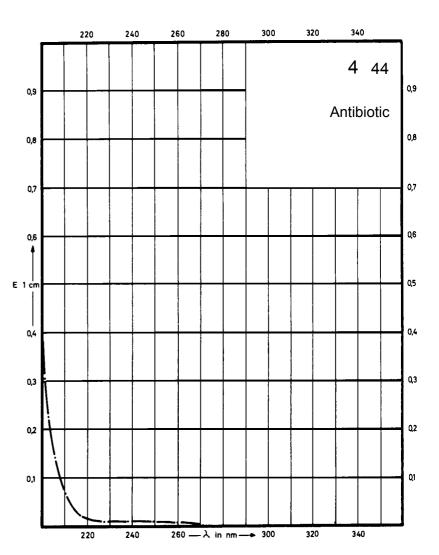

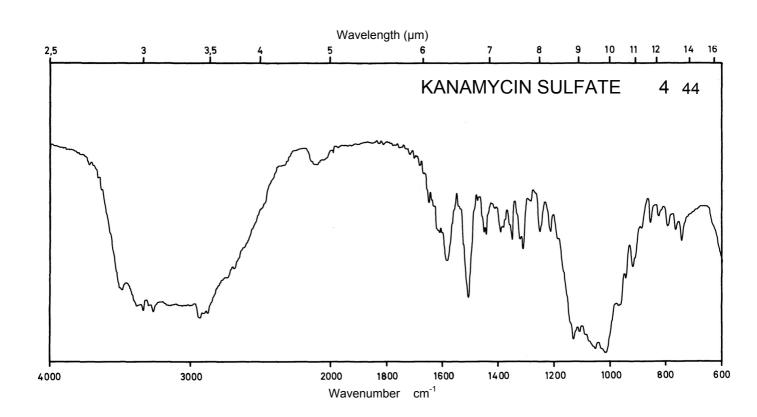



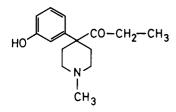



ε

5 mg / 100 ml Concentration


247.3

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of absorption 280 nm 94 2320

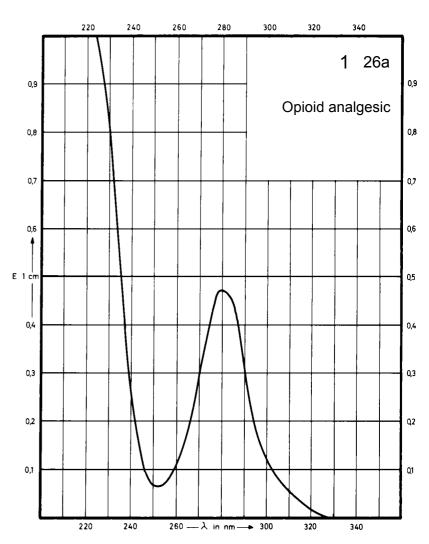

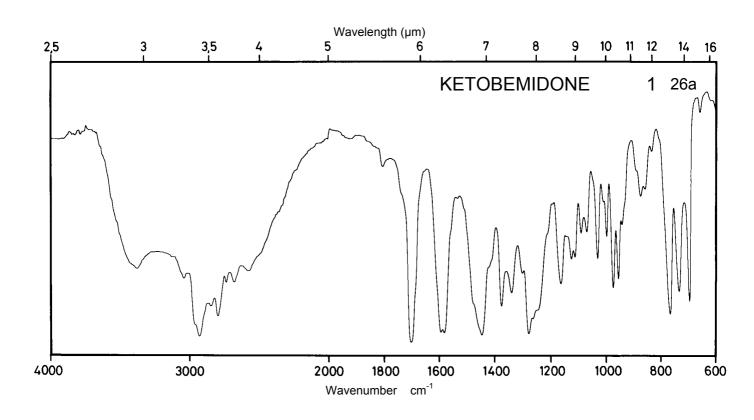



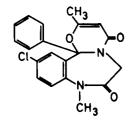


 M_{r}


368.8

Concentration


0.85 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm	239 nm	Decom- position observed	239 nm
E ^{1%} 1cm	490	514		524
ε	18070	18970		19320

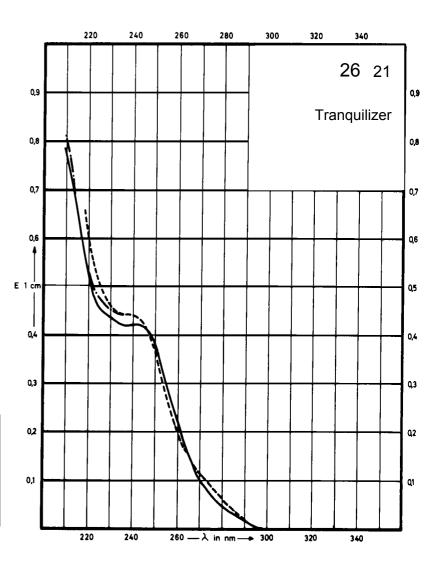

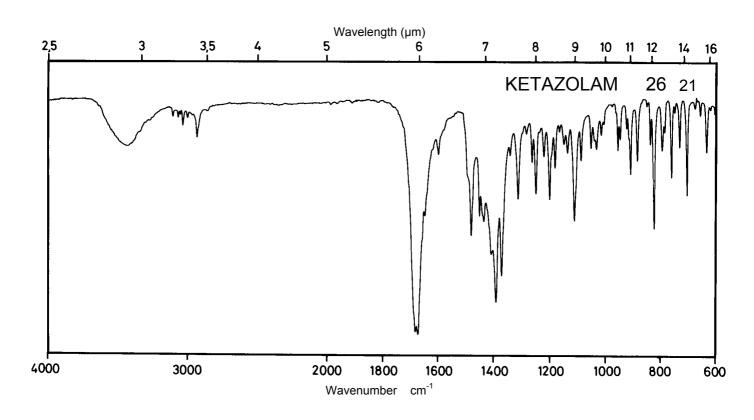




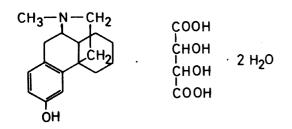



443.5


Concentration 2 mg / 100 ml 6 mg / 100 ml 8 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm		279 nm	299 nm 239 nm
E ^{1%} 1cm	52		45	68 195
ε	2290		2010	3030 8650

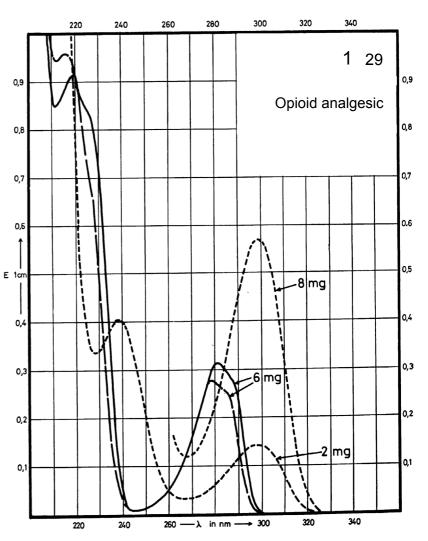


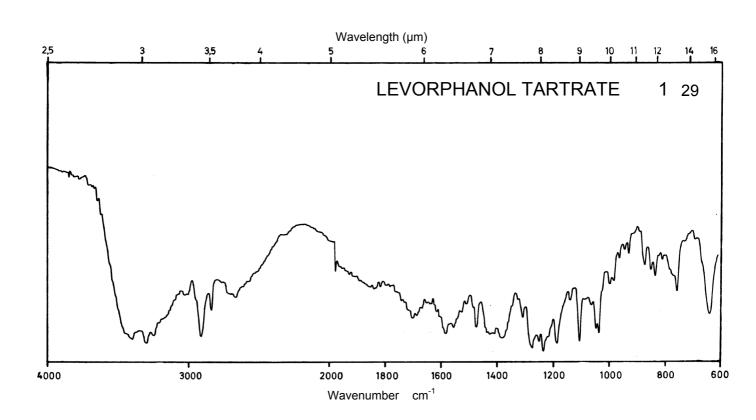

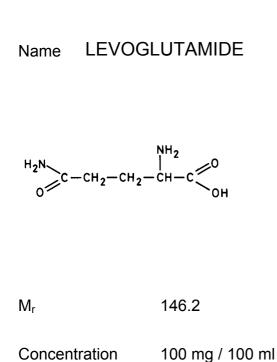


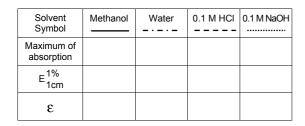

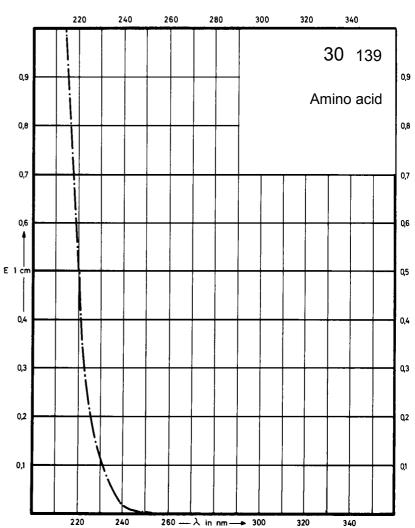

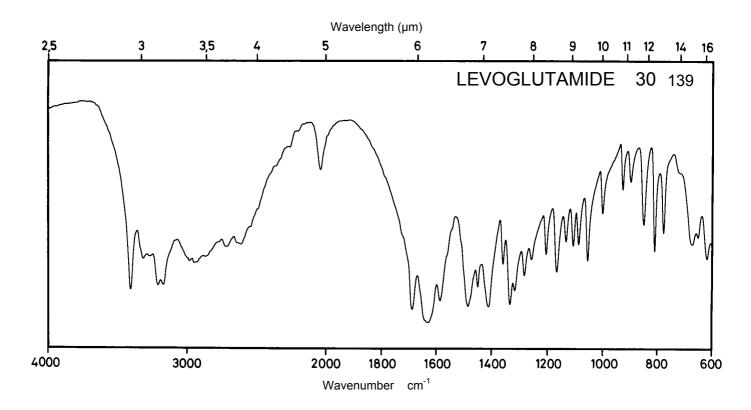





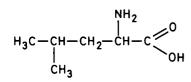



131.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

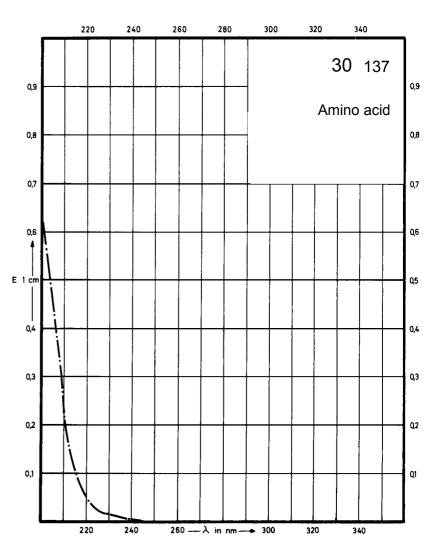


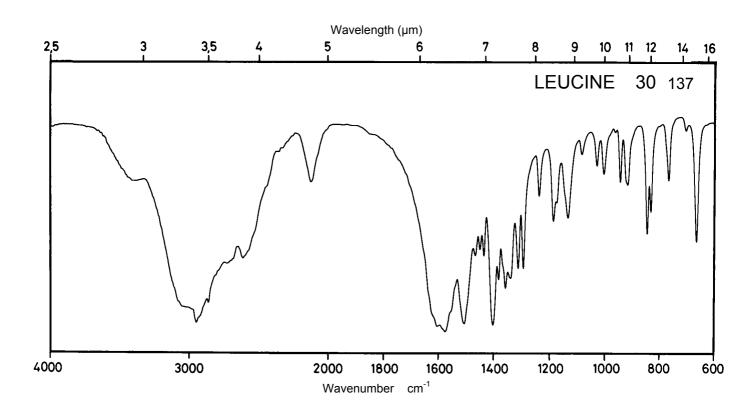

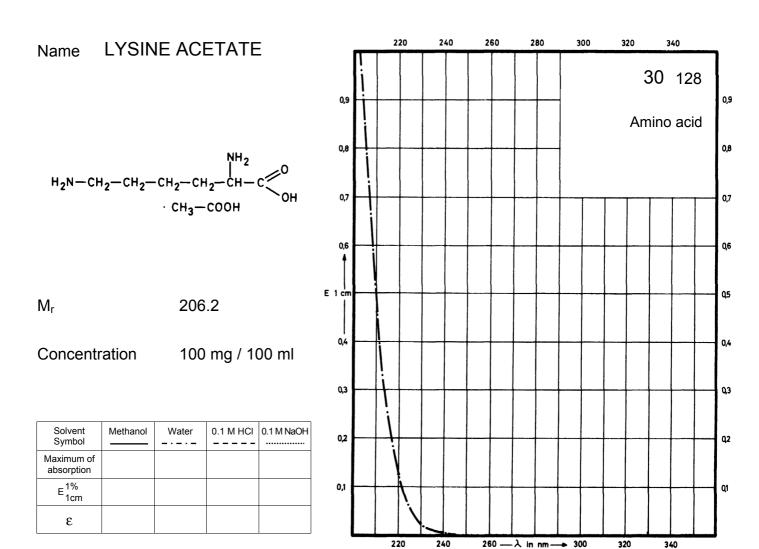

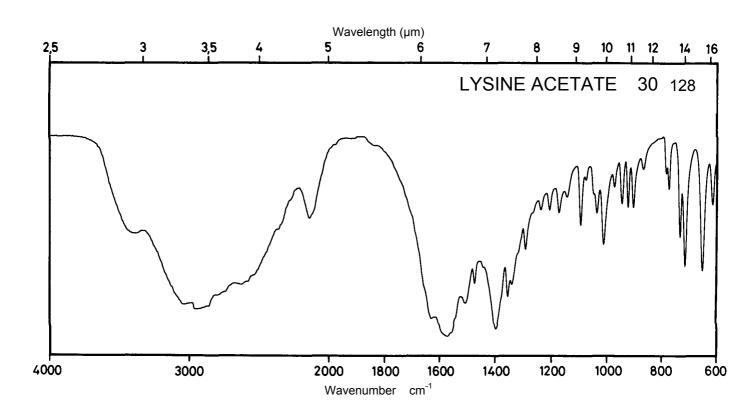





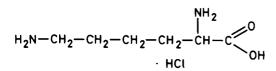


182.7


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
ε				

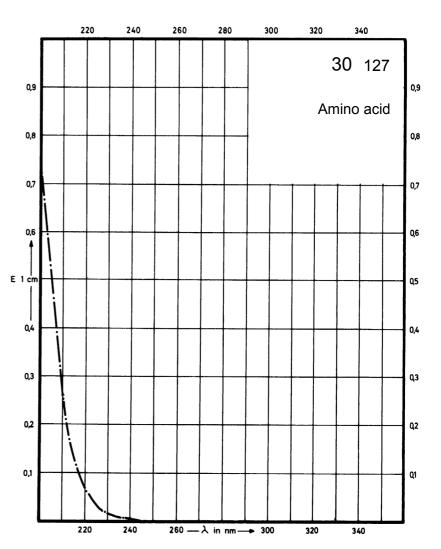


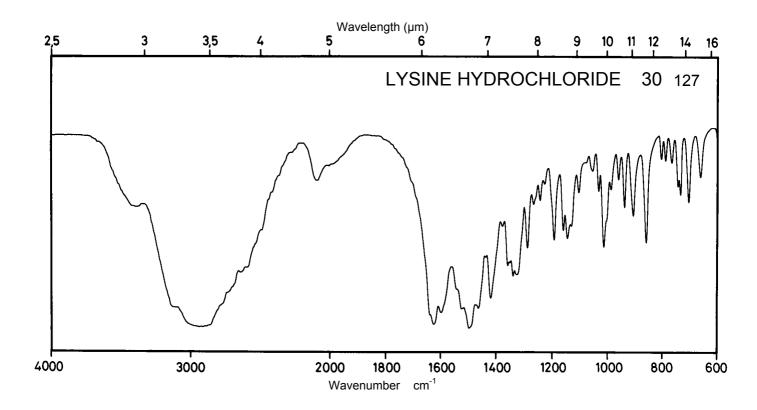

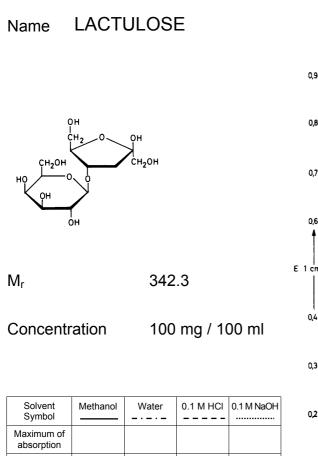

3

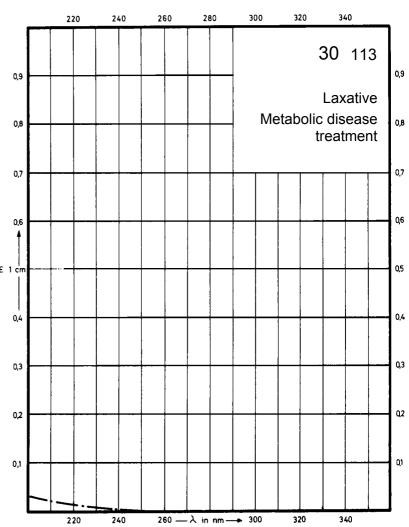

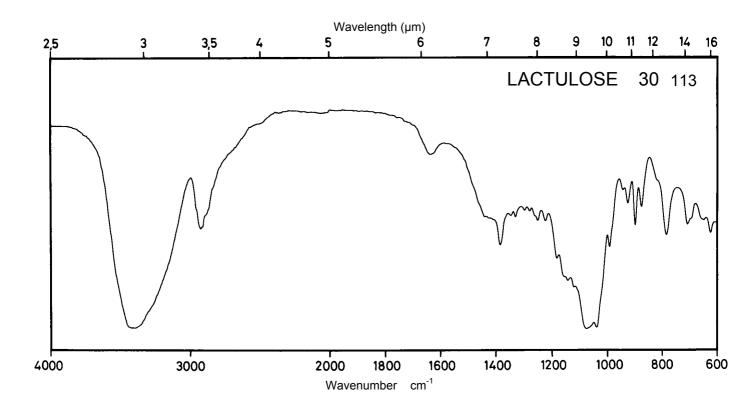








 M_r 290.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

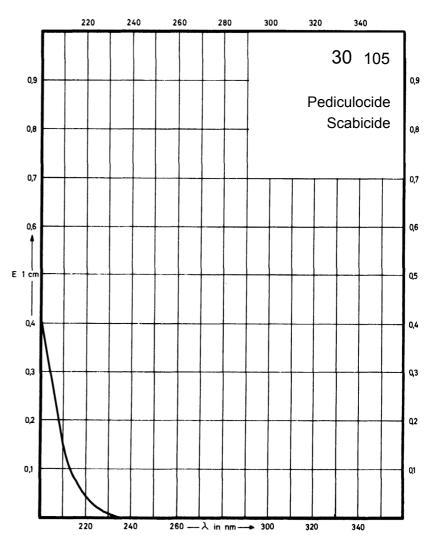

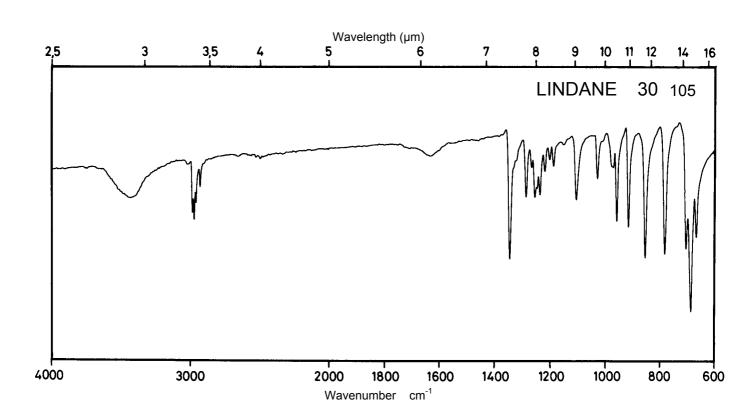





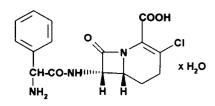



367.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm	262 nm	269 nm	Decom- position observed
E ^{1%} 1cm	309	320	290	
ε	11400	11800	10700	

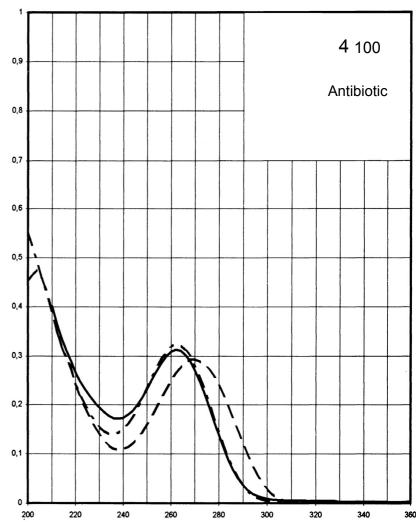

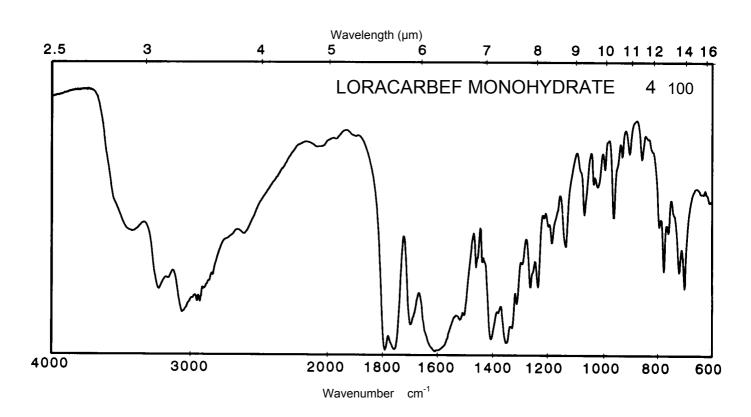

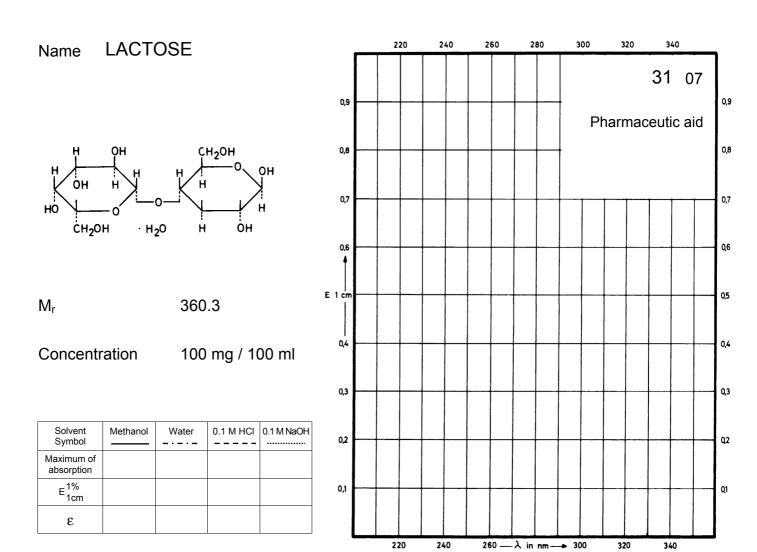

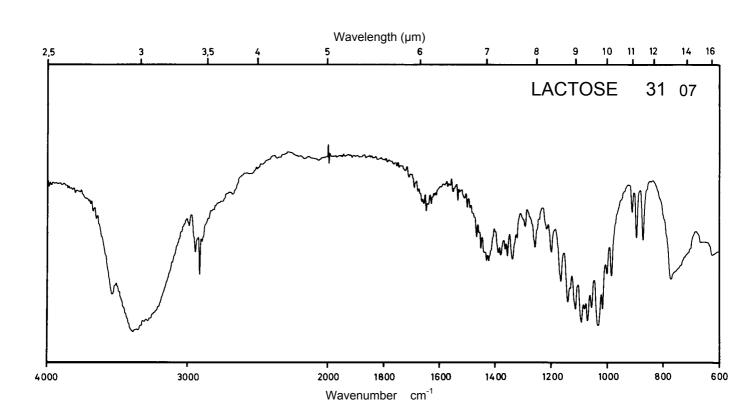



Name LISURIDE MALEATE


 M_{r}

454.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	313 nm	310 nm	312 nm	309 nm
	240 nm	238 nm	238 nm	237 nm
E ^{1%} 1cm	190	182	182	199
	487	504	489	557
ε	8640	8270	8270	9040
	22130	22900	22230	25320

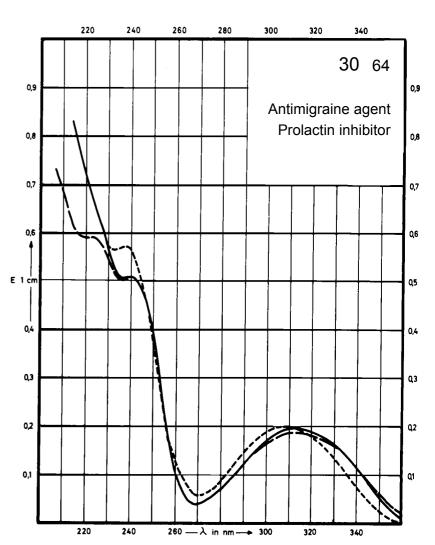

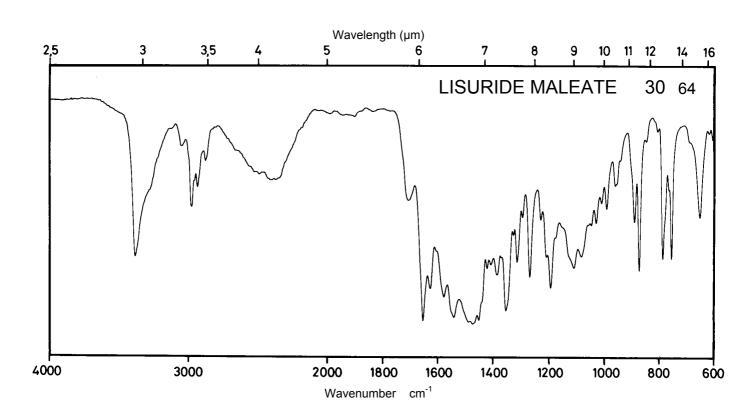

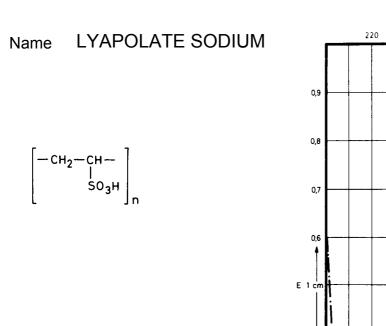

Concentration


Methanol

Water

Solvent

Symbol Maximum of absorption

3

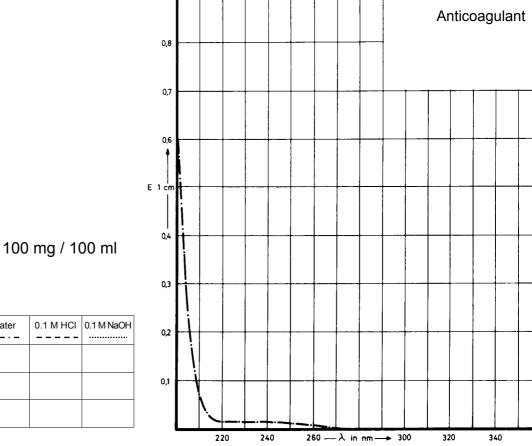

240

280

260

300

320

340

30 32

0,8

0,7


0,6

0,5


0,4

0,3

Q2

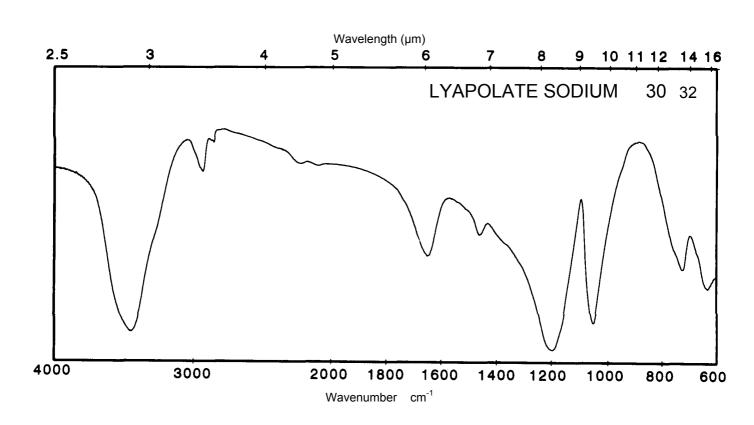



Name LEVALLORPHAN TARTRATE

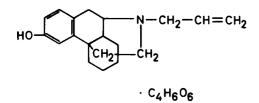


 M_r 433.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		278 nm	298 nm
E ^{1%} 1cm	54		47	69
ε	2340		2040	2990

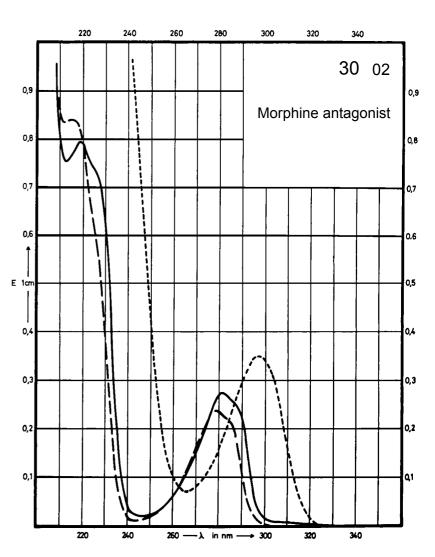

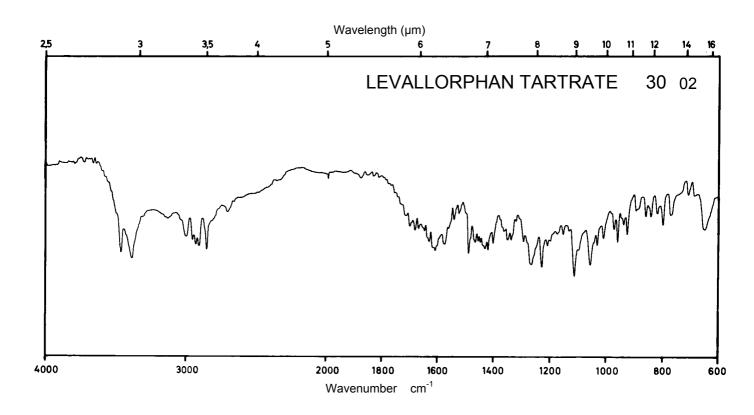





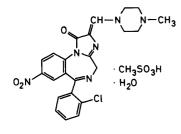


579.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	330 nm	329 nm	Decom- position observed	Decom- position observed
E ^{1%} 1cm	700	675		
ε	40600	39100		

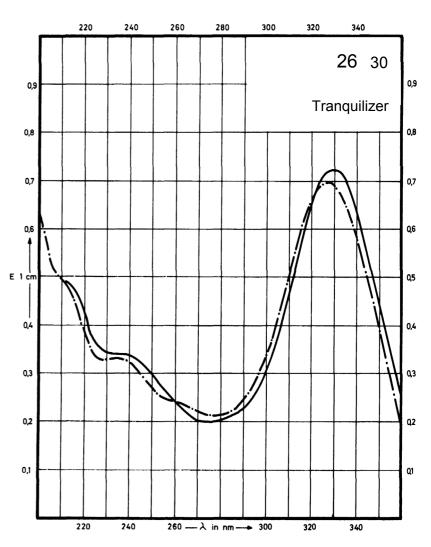

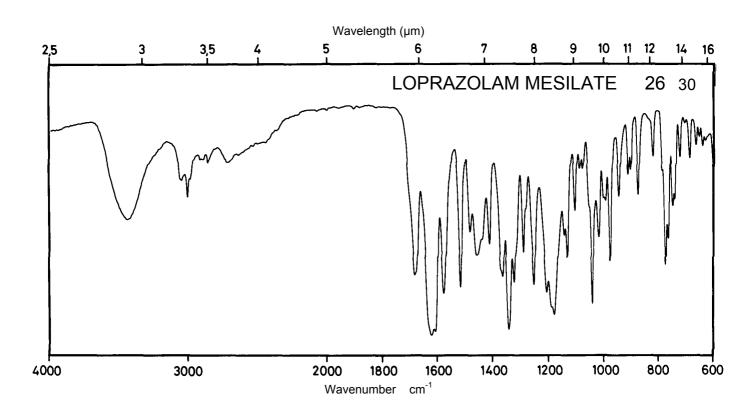




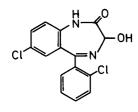


321.2


Concentration 0.5 mg / 100 ml 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	319 nm 229 nm		231 nm	346 nm
E ^{1%} 1cm	64 1160		1110	82
3	2050 37260		35550	2650

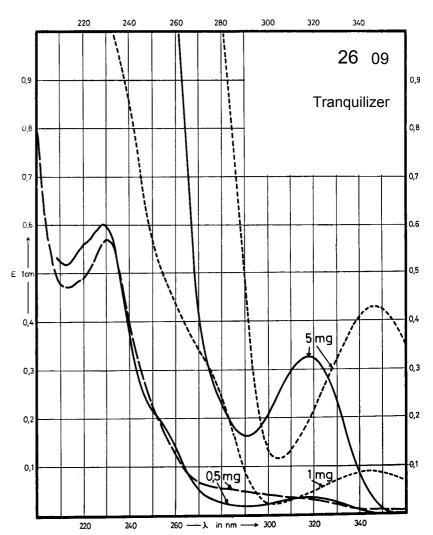

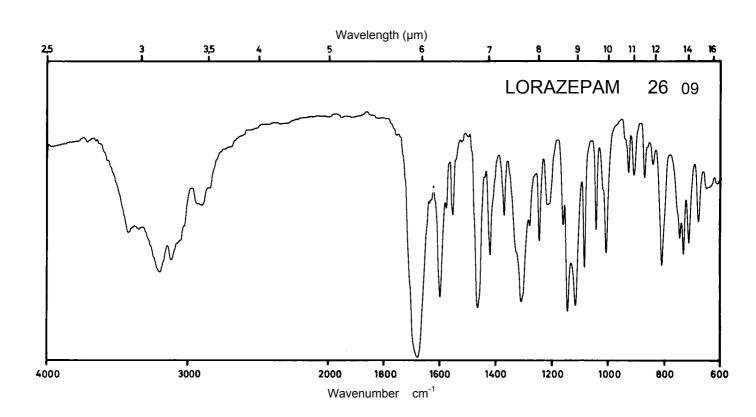





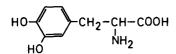

197.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			279 nm	
E ^{1%} 1cm			141	
ε			2780	

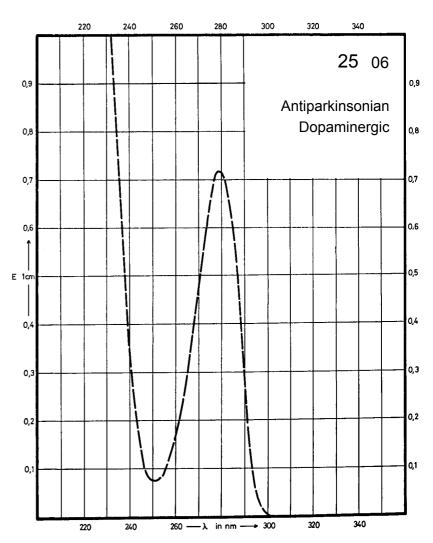

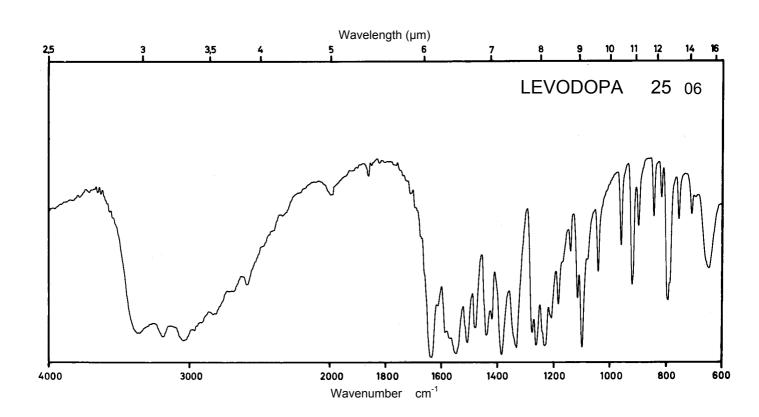





Glp-His-Trp-Ser-Tyr-D -Leu-Leu-Arg-Pro-NHC $_2$ H $_5$ x CH $_3$ COOH


 M_{r} 1269.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm	275 nm	274 nm	277 nm
E ^{1%} 1cm	60	49	49	52
3	7610	6260	6210	6550

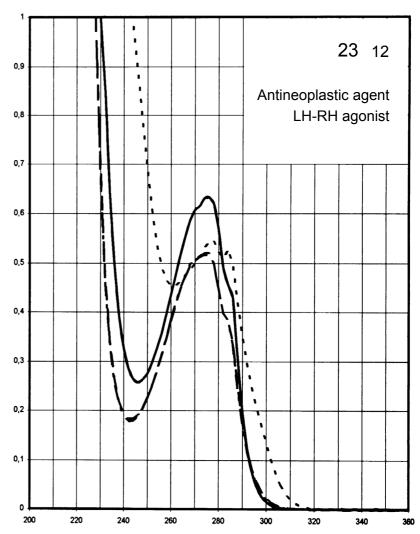

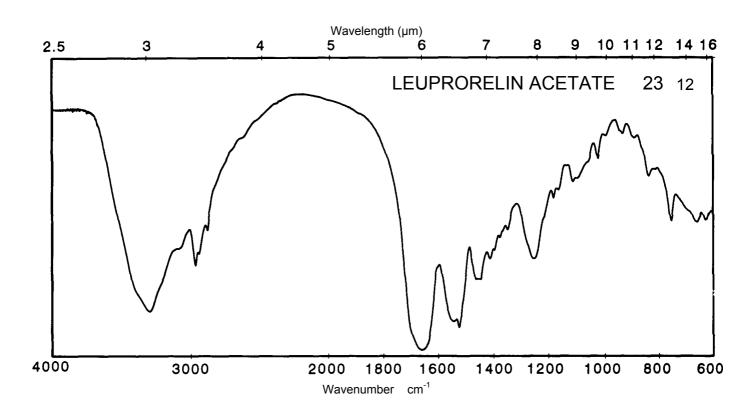



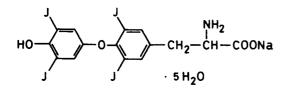



889.0

Concentration


8.4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				325 nm
E ^{1%} 1cm				68
ε				6000

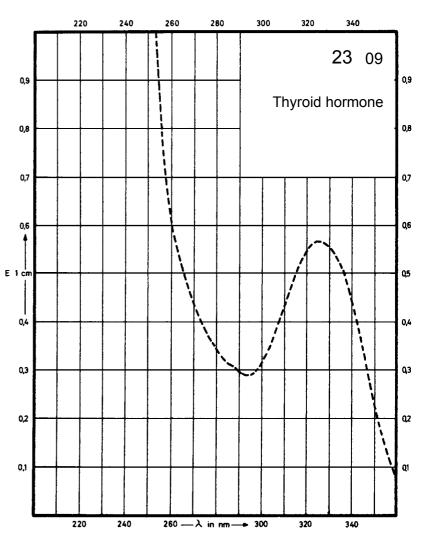

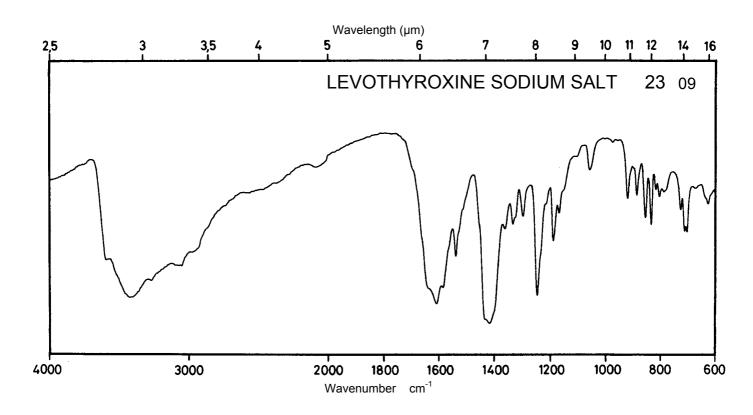




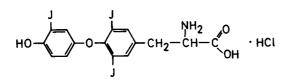



687.4

Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	298 nm 226 nm		294 nm 224 nm	319 nm
E ^{1%} 1cm	59 600		56 557	64
3	4060 41240		3850 38290	4400

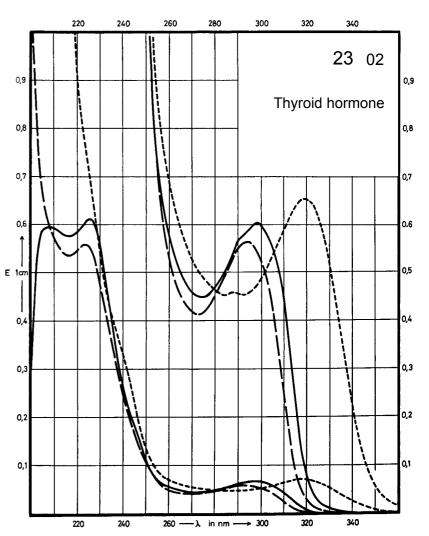


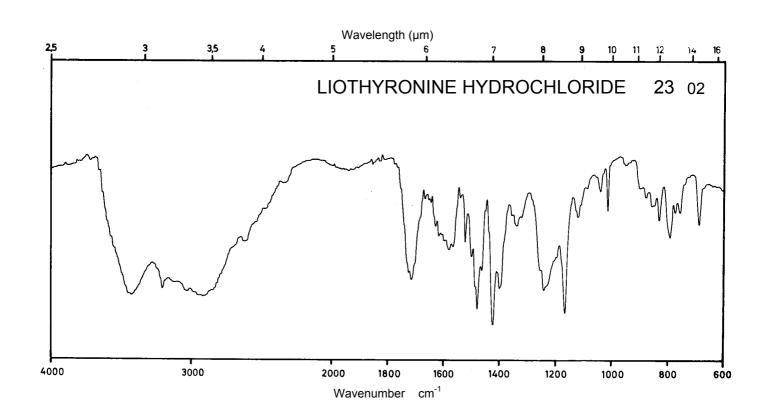


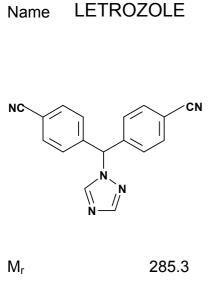


Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm	281 nm	282 nm	281 nm
	273 nm	273 nm	274 nm	273 nm
	240 nm	241 nm	240 nm	241 nm
E ^{1%} 1cm	55 59 1177	63 69 1192	66 72 1195	62 69 1196
3	1560	1790	1870	1780
	1680	1960	2060	1970
	33600	34000	34100	34100

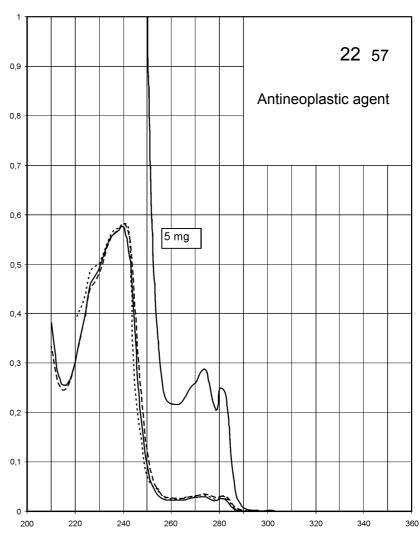

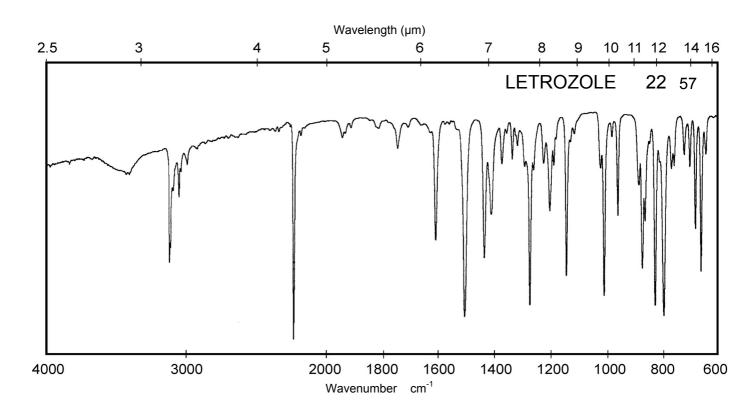




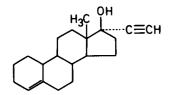


 M_r 284.4


Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

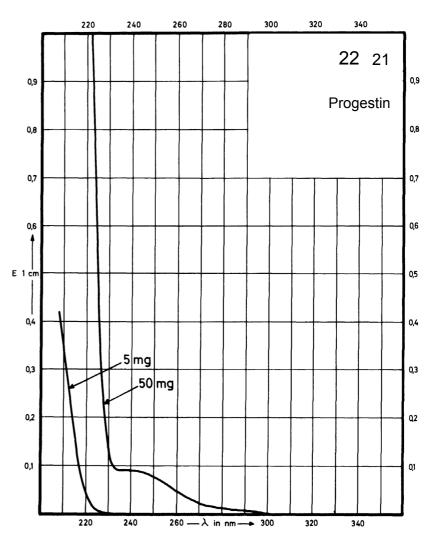

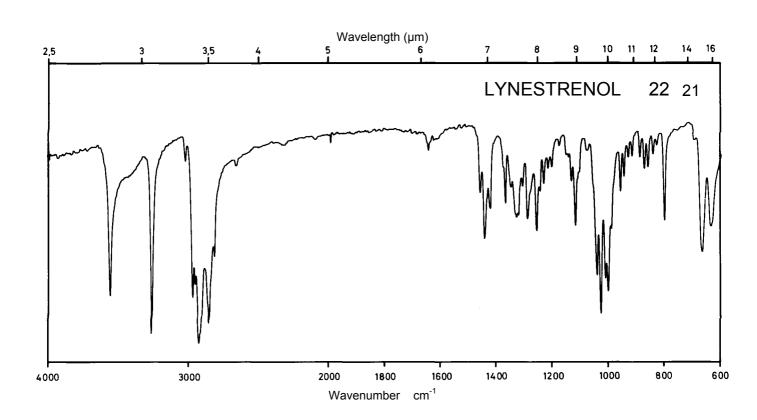





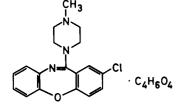


 M_r 445.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	320 nm 297 nm 249 nm		293 nm 252 nm	298 nm 249 nm
E ^{1%} 1cm	149 153 284		218 251	188 273
ε	6650 6840 12660		9700 11200	8370 12170

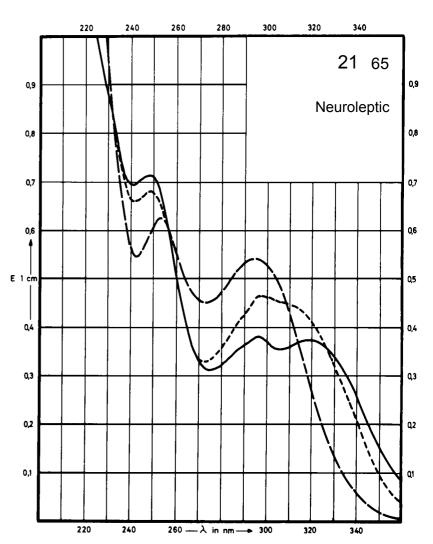

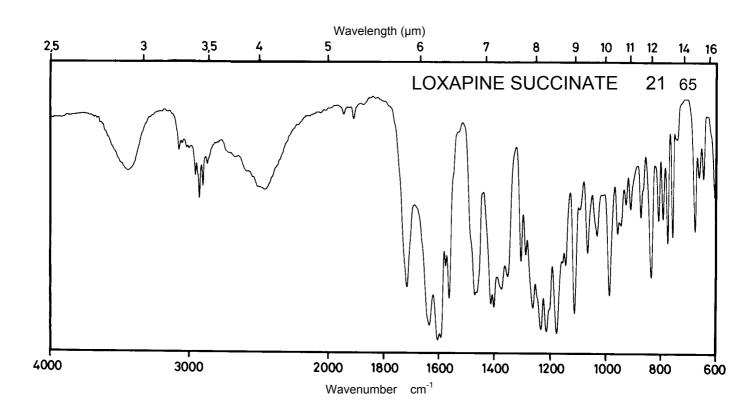




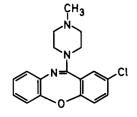


 M_r


327.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	320 nm 298 nm 251 nm		293 nm 252 nm	298 nm 249 nm
E ^{1%} 1cm	210 214 390		299 346	259 374
ε	6870 7010 12800		9810 11330	8490 12270

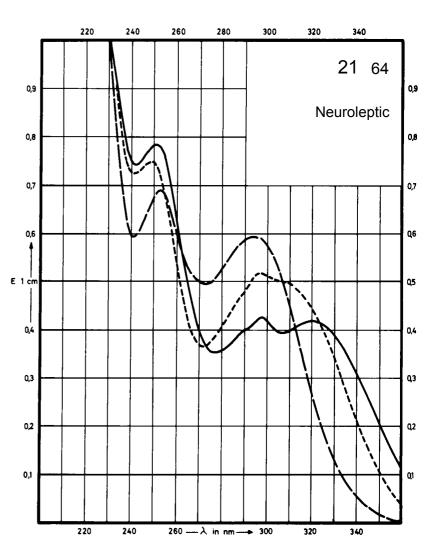

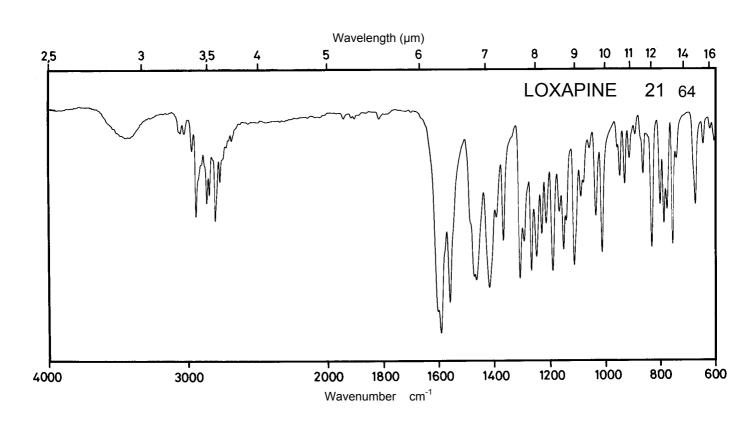





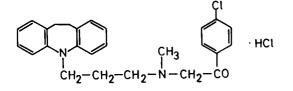

 M_r 455.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	255 nm	262 nm	262 nm	
E ^{1%} 1cm	545	441	521	
ε	24800	20060	23750	

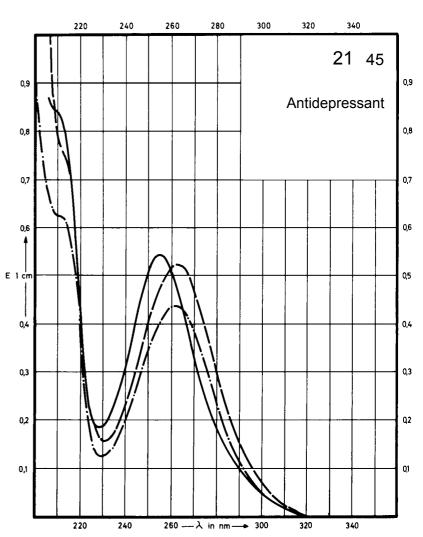

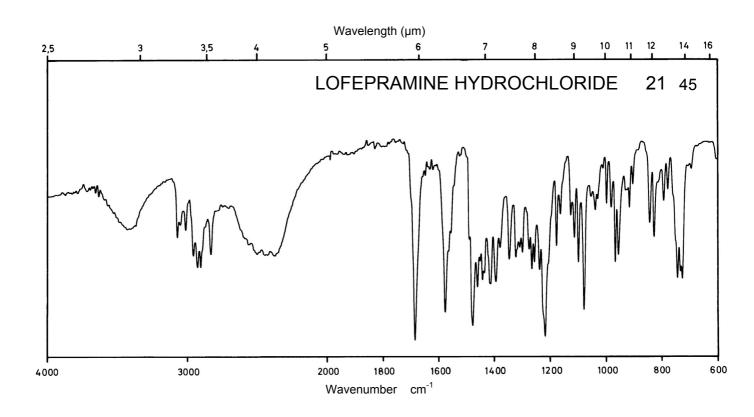


Name LEVOMEPROMAZINE MALEATE


 M_{r}

444.6


Concentration 1 mg / 100 ml 5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	304 nm 253 nm		302 nm 250 nm	
E ^{1%} 1cm	100 595		93 595	
ε	4450 26450		4130 26450	

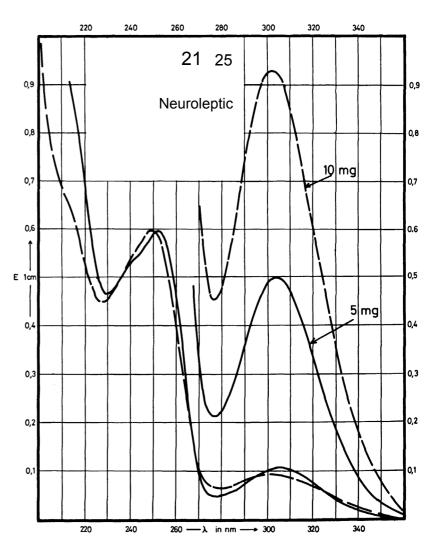

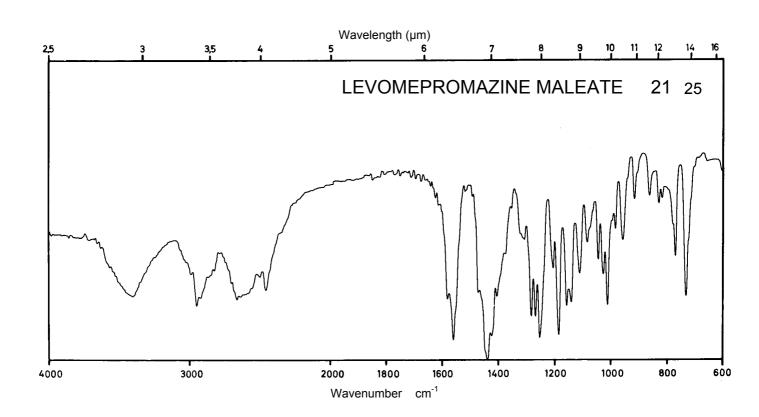



Name LEVOBUNOLOL HYDROCHLORIDE

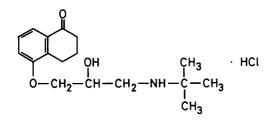

 M_r 327.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm 254 nm		315 nm 257 nm	316 nm 258 nm
E ^{1%} 1cm	70 266		68 283	67 277
ε	2300 8700		2200 9300	2200 9100

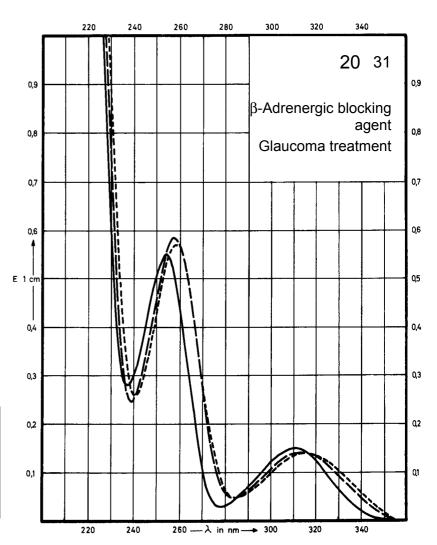

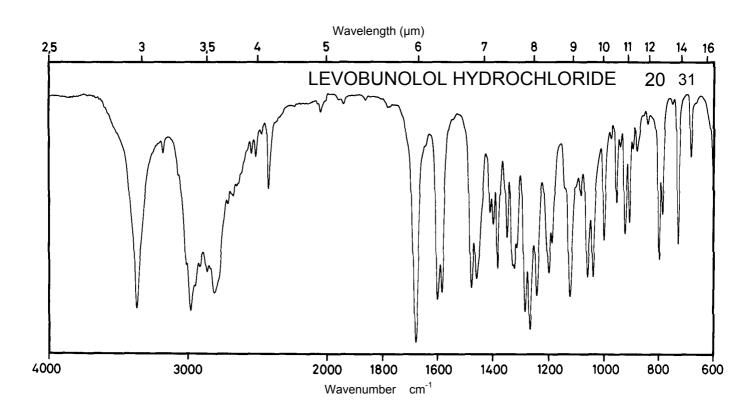





 M_{r} 364.9


Concentration 2.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	305 nm		303 nm	333 nm 246 nm
E ^{1%} 1cm	96		86	145 241
3	3490		3140	5290 8800

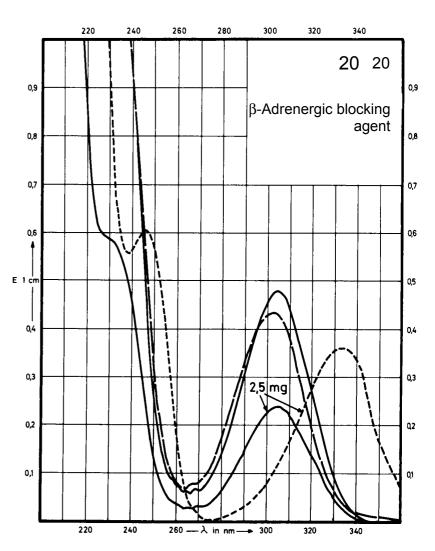

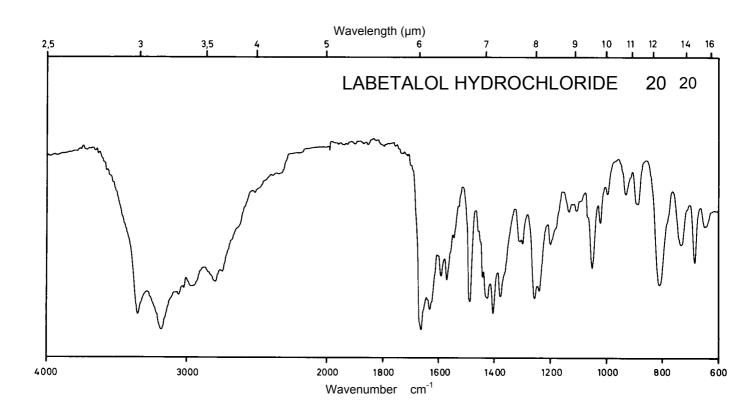





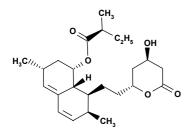

 M_r 404.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	246 nm	247 nm	247 nm	248 nm
	238 nm	239 nm	239 nm	239 nm
	230 nm	232 nm	232 nm	231 nm
E ^{1%} 1cm	423 631 545	369 555 498	359 541 485	369 558 518
3	17100	14900	14500	14900
	25500	22500	21900	22600
	22100	20100	19600	21000

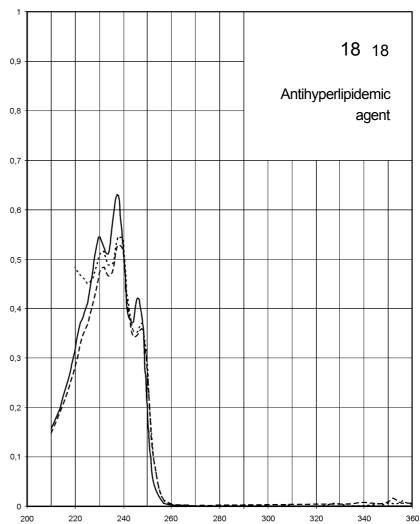

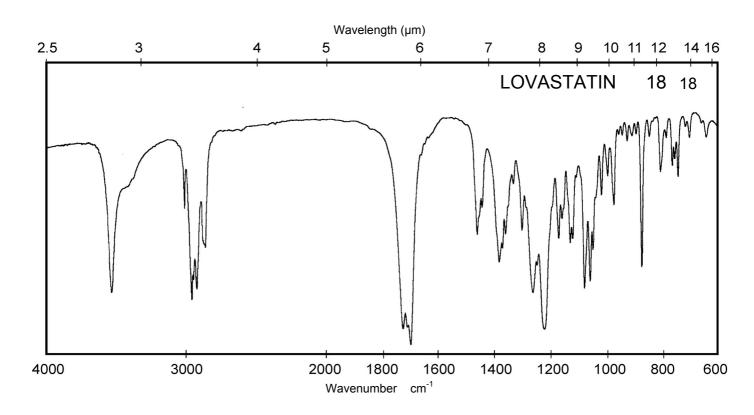



Name LIDOCAINE HYDROCHLORIDE MONOHYDRATE

$$\begin{bmatrix} CH_{3} & & & \\ -NH-CO-CH_{2}-N & C_{2}H_{5} \\ CH_{3} & & CI-H_{2}O \end{bmatrix}$$


 M_r 288.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 263 nm	271 nm 263 nm	271 nm 263 nm	
E ^{1%} 1cm	10.3 13.4	12.8 15.7	12.8 15.7	
ε	297 387	370 453	370 453	

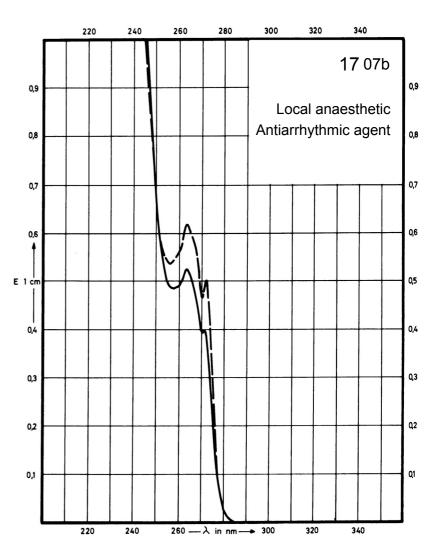

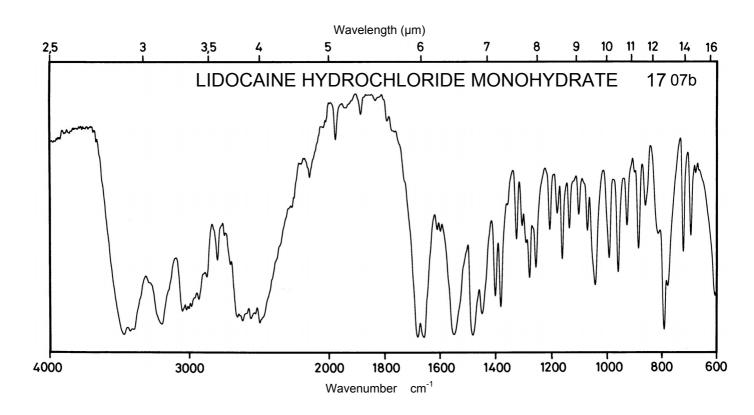




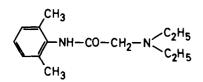


234.3


30 mg / 100 ml Concentration

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 263 nm		271 nm 263 nm	
E ^{1%} 1cm	10.3 14.8		15.9 19.3	
ε	240 347		373 453	

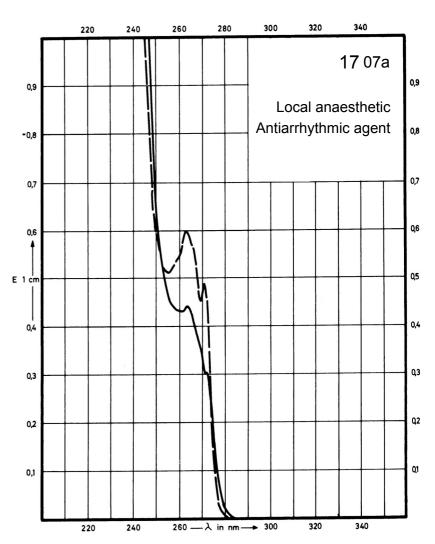

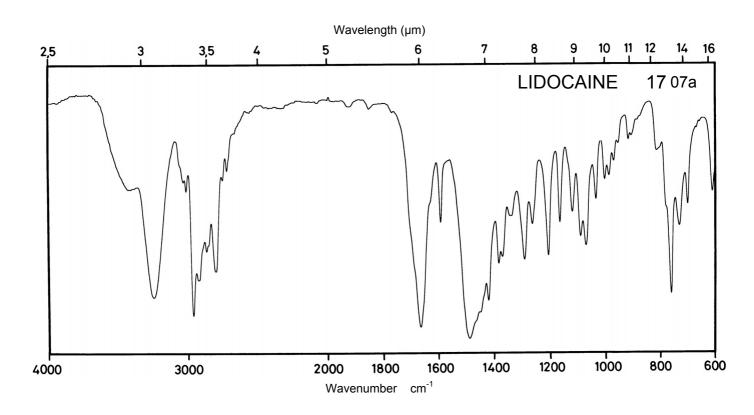



Name LIDOCAINE HYDROCHLORIDE

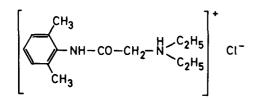


 M_r 270.8

			1	
Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
-,				
Maximum of absorption	271 nm 263 nm	271 nm 263 nm	271 nm 263 nm	
E ^{1%} 1cm	10.9 14.2	14.1 16.8	14.1 16.8	
ε	295 385	382 455	382 455	

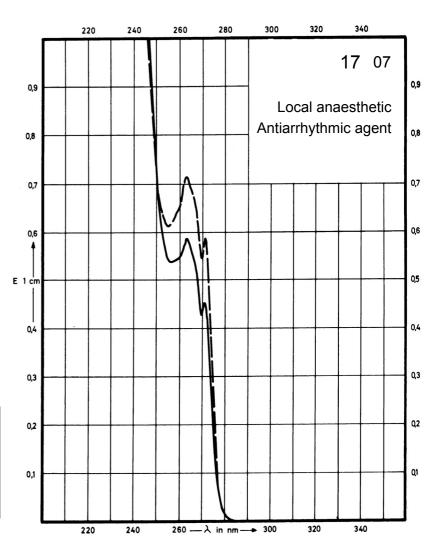

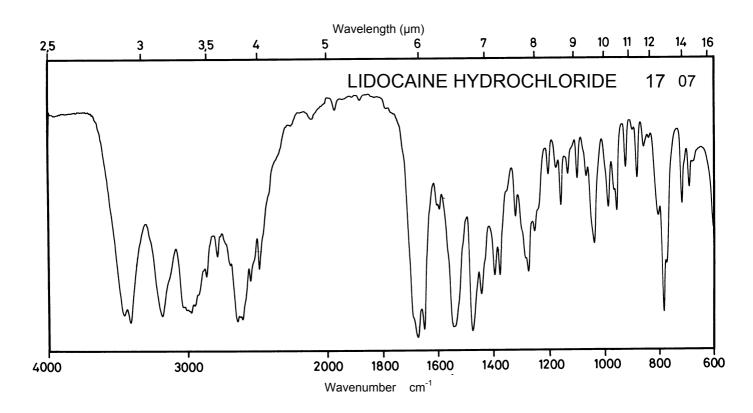





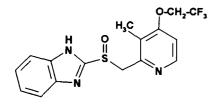


369.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	Decom- position observed	Decom- position observed	292 nm
E ^{1%} 1cm				466
ε				17200

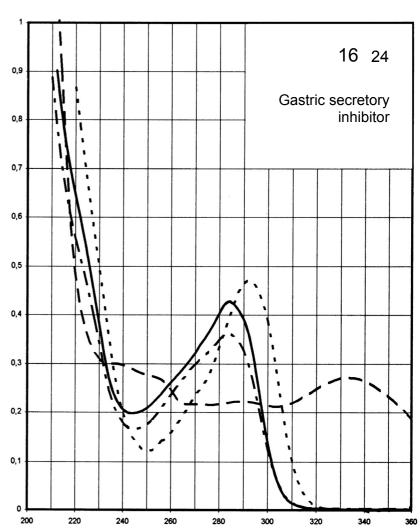

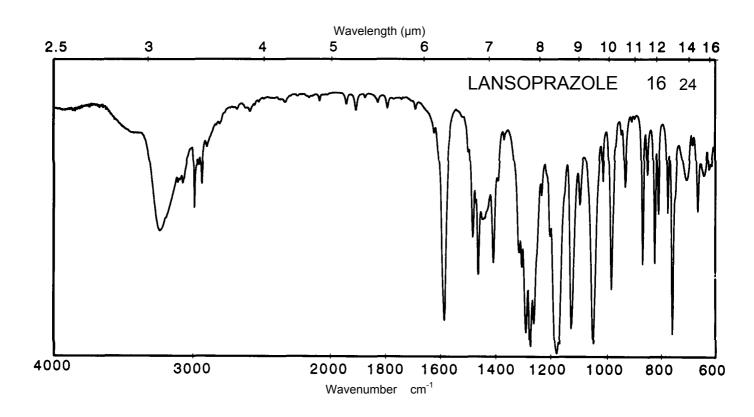





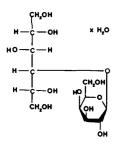


 M_r 362.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	227 nm			
E ^{1%} 1cm	4.7			
ε	170			

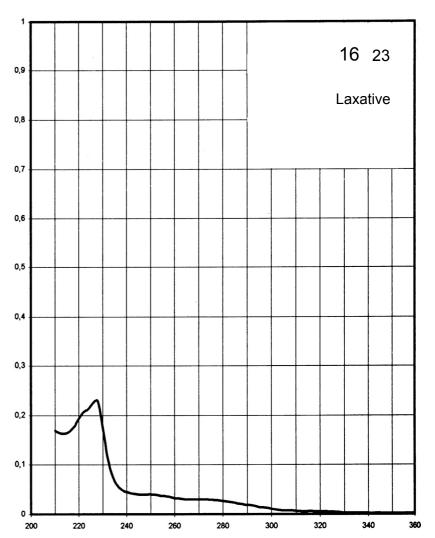

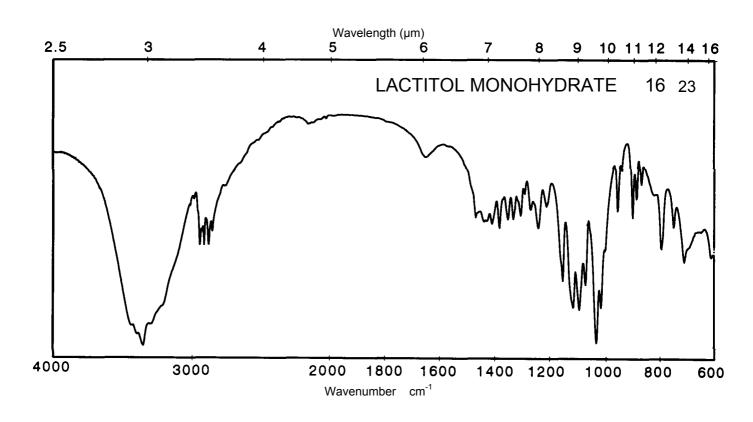





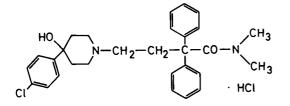


513.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm 260 nm			
E ^{1%} 1cm	11.4 12.8			
ε	587 655			

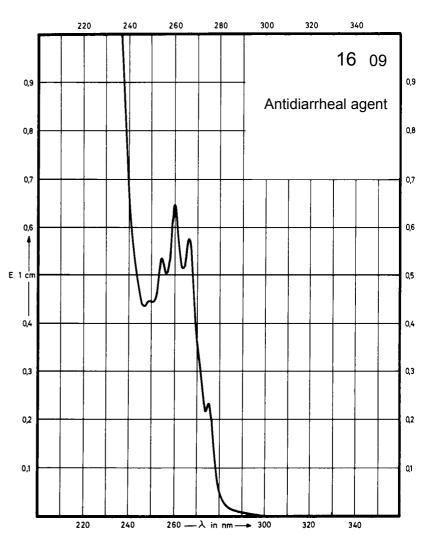

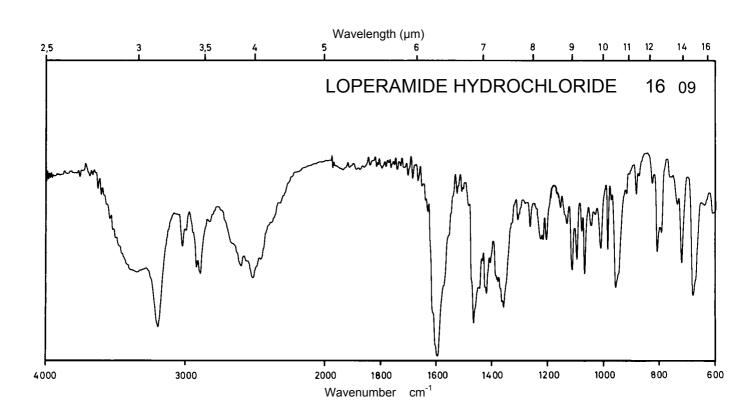




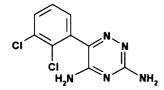


 M_r


256.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	Decom- position observed	267 nm	307 nm
E ^{1%} 1cm			277	285
3			7100	7300

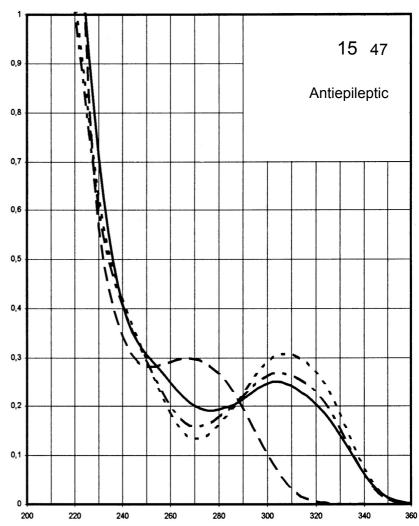

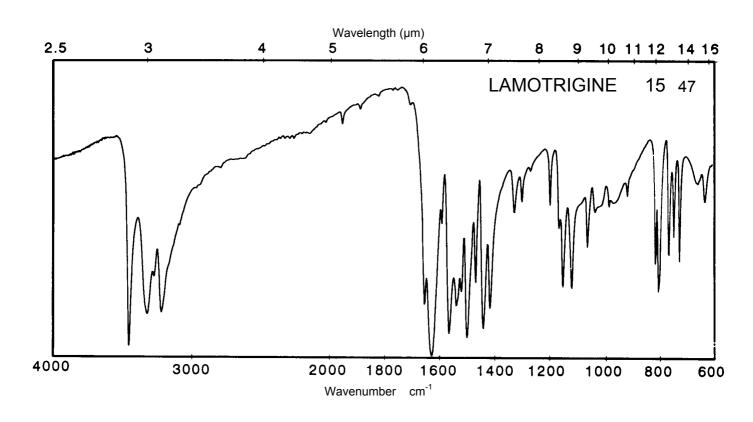




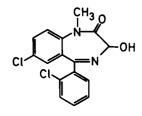


 M_r 335.2


Concentration 0.5 mg / 100 ml 10 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 316 nm 230 nm 314 nm 231 nm 311 nm 231 nm 313 nm 229 nm absorption E^{1%} 1cm 55 1000 54 1070 59 1030 54 1020 1840 33520 1820 35870 1980 34520 1820 34130 ε

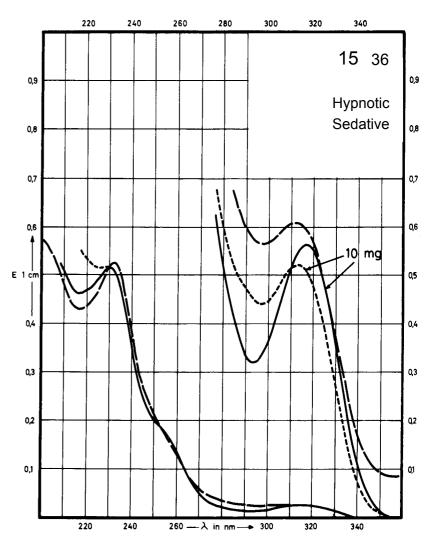

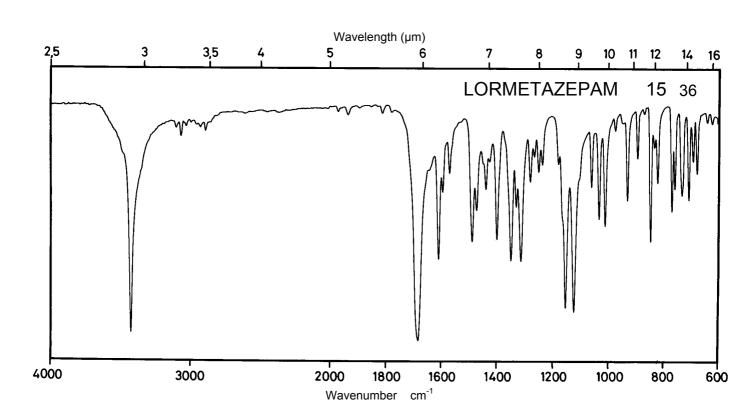




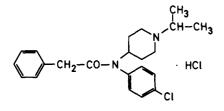


 M_r 407.4


Concentration 2 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 259 nm			
E ^{1%} 1cm	11 13			
ε	448 525			

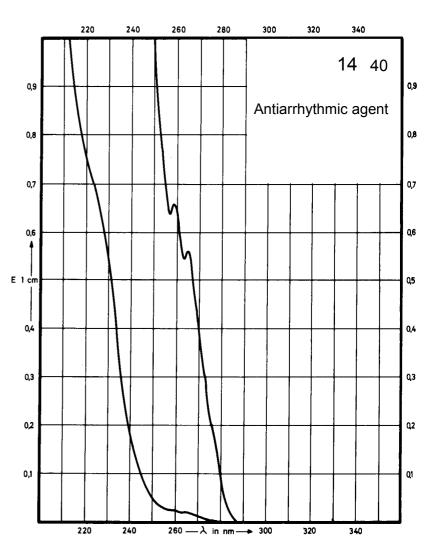

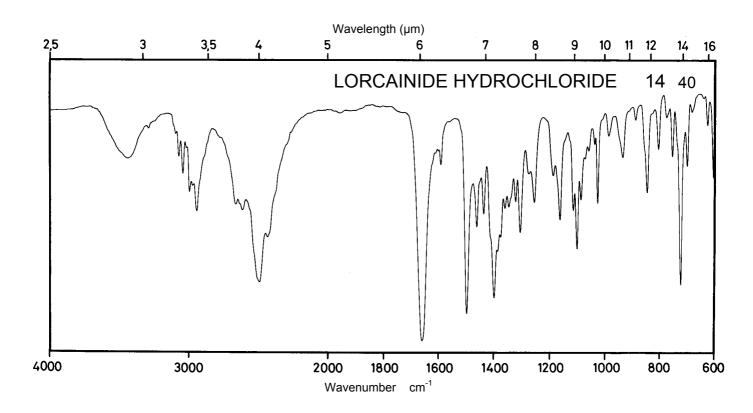





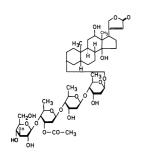


 M_r 985

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	219 nm			
E ^{1%} 1cm	150			
3	14780			

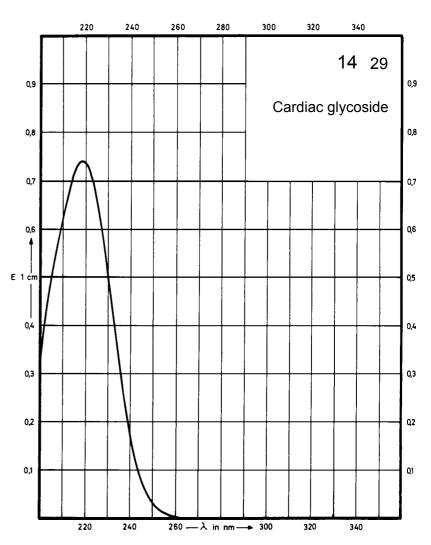

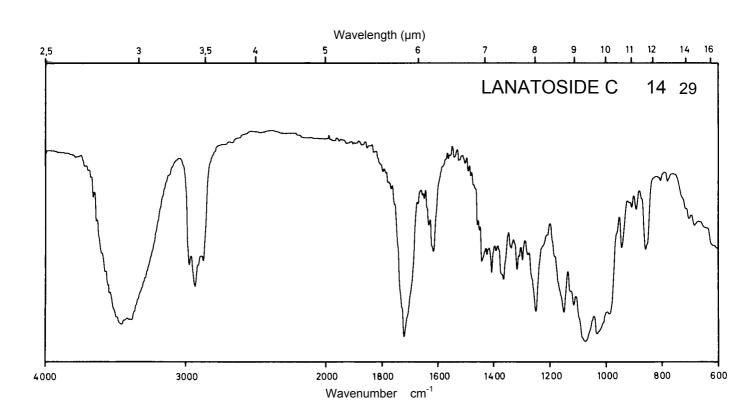





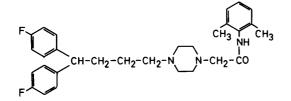


491.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 266 nm			272 nm 265 nm
E ^{1%} 1cm	39 43			36 42
3	1906 2120			1770 2060

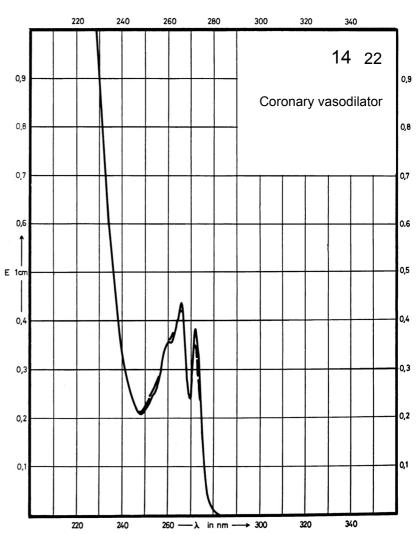

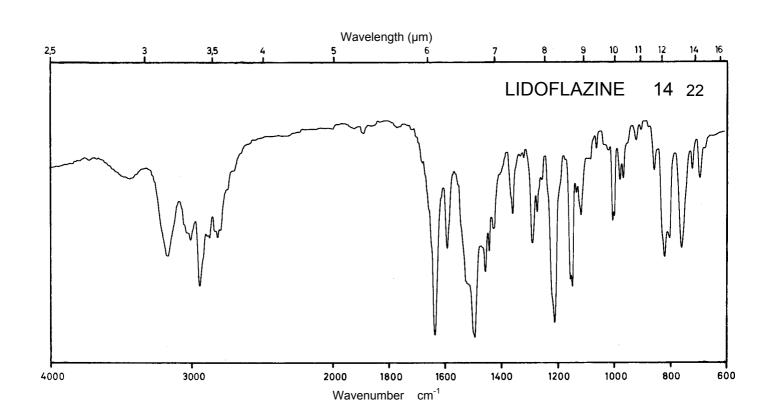





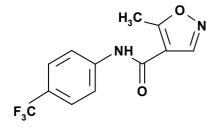


 M_r 270.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm	259 nm	259 nm	Decom- position observed
E ^{1%} 1cm	732	610	610	
ε	19800	16500	16500	

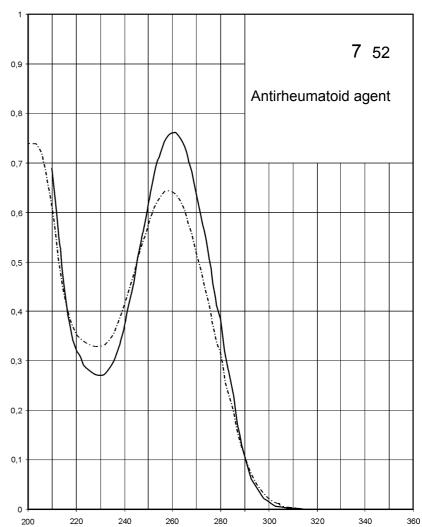

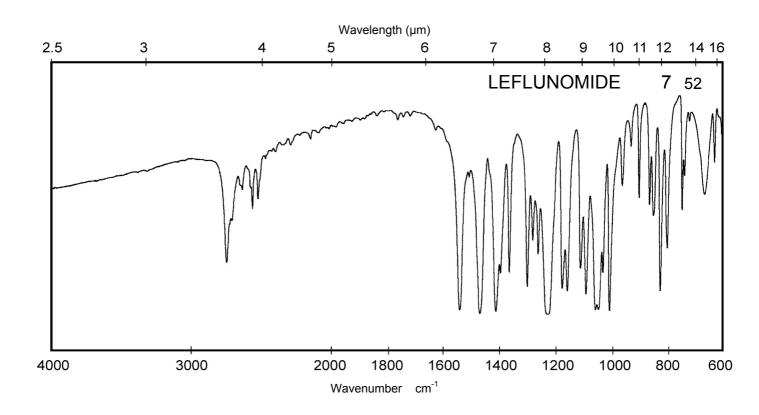


Name LONAZOLAC CALCIUM SALT

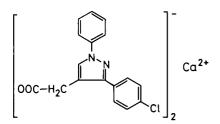



 M_{r}


663.6

Solvent Symbol	Methanol	Water _ · - · -	0.1 M HCI	0.1 M NaOH
Maximum of absorption			274 nm	279 nm
E ^{1%} 1cm			652	700
ε			43270	46450

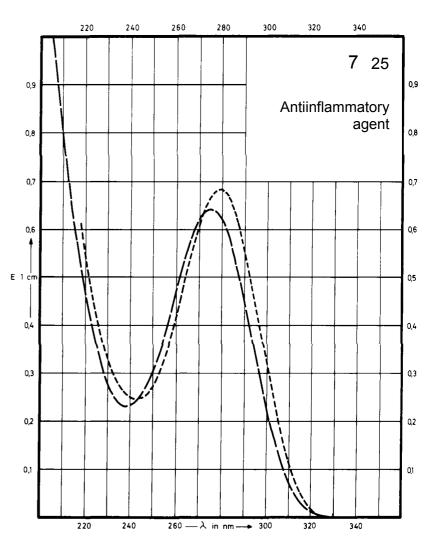


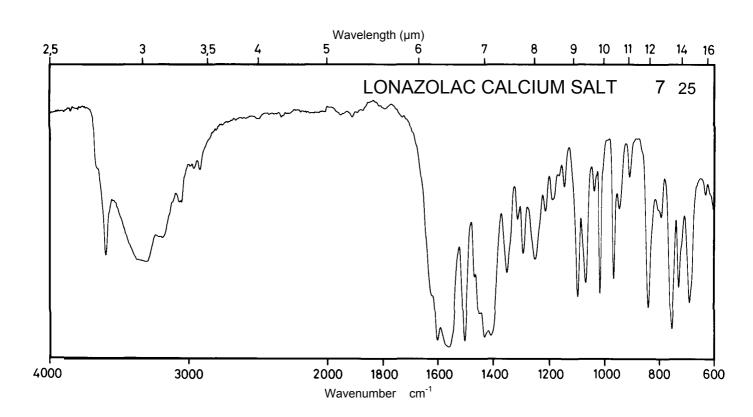

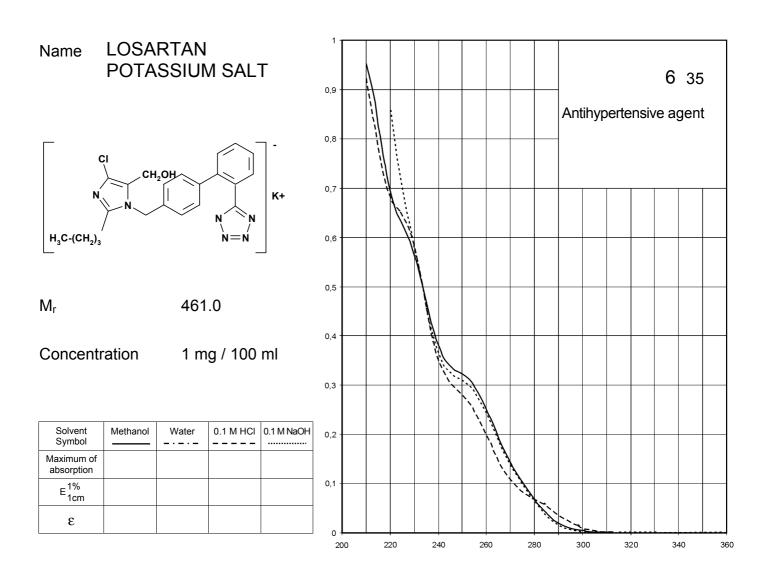

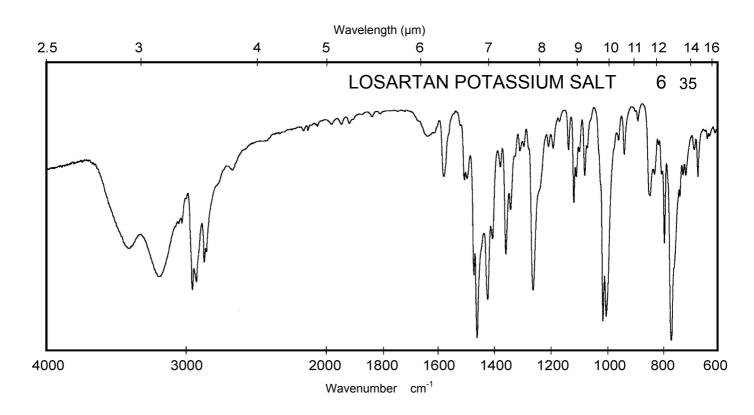



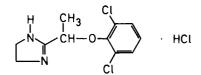



 M_r


Concentration 1 mg / 100 ml 50 mg / 100 ml


295.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
-,				
Maximum of absorption	270 nm			
E ^{1%} 1cm	10.0			
3	2950			

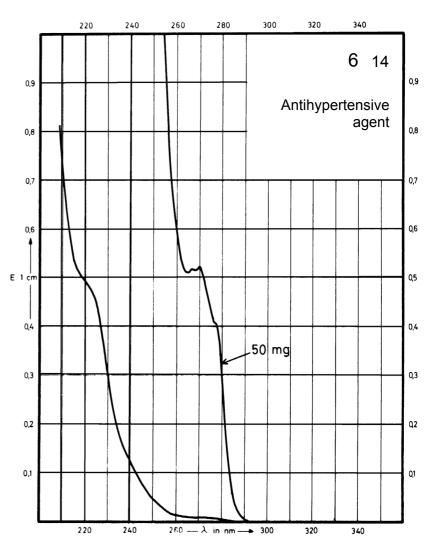


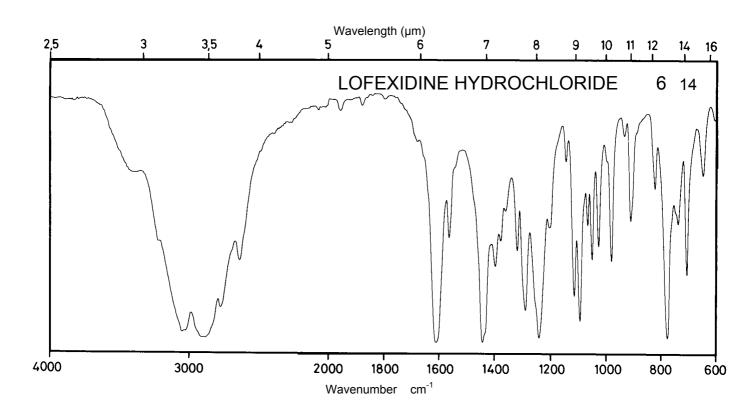

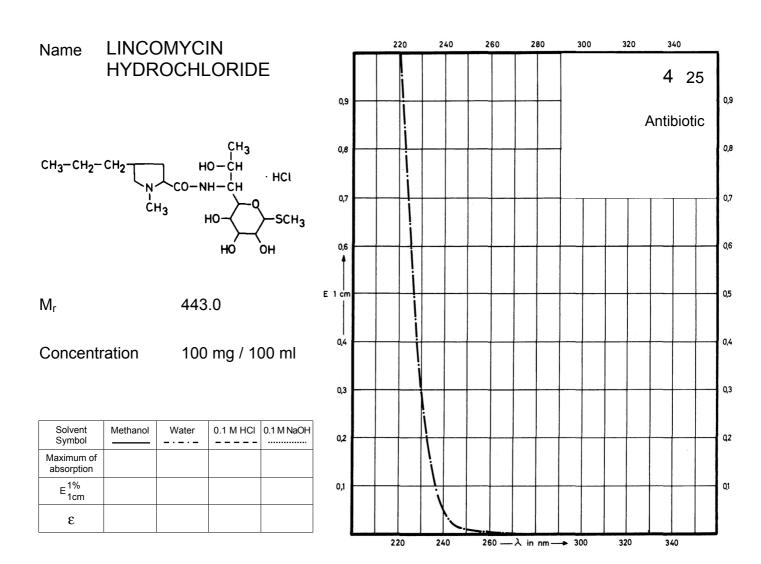

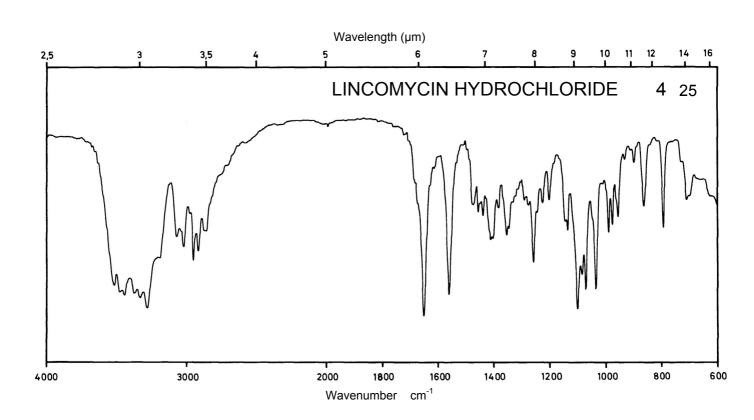



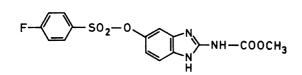



365.3

Concentration


0.84 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	290 nm 285 nm 245 nm		285 nm 279 nm 221 nm	
E ^{1%} 1cm	385 390 305		480 435 985	
ε	14060 14250 11150		17500 15900 36000	

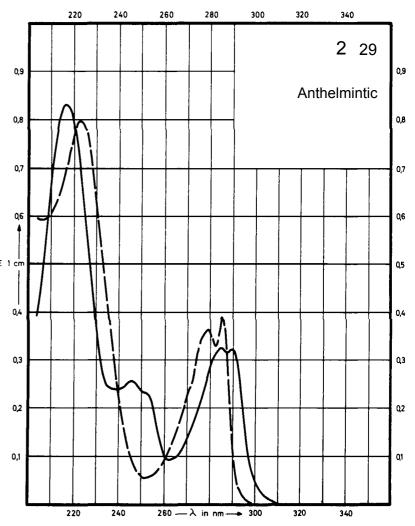

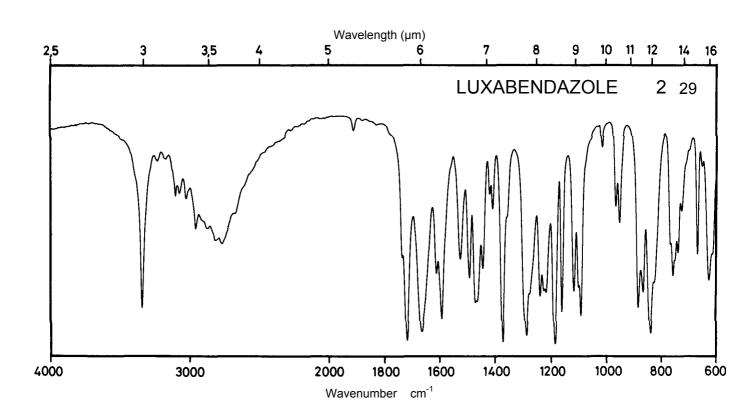



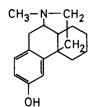



257.4


Concentration


8.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm			
E ^{1%} 1cm	90			
ε	2300			

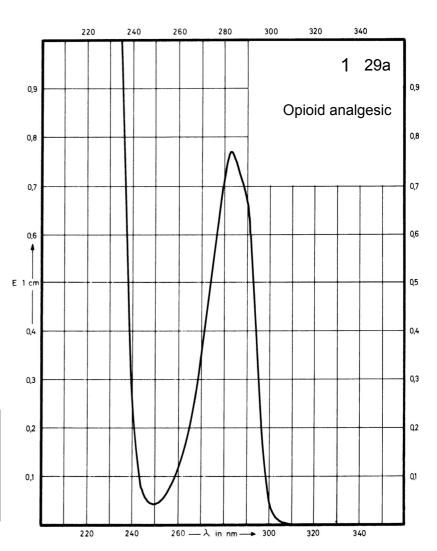


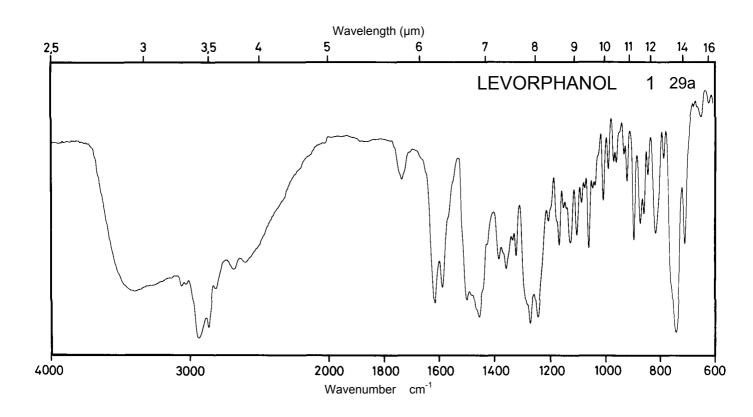

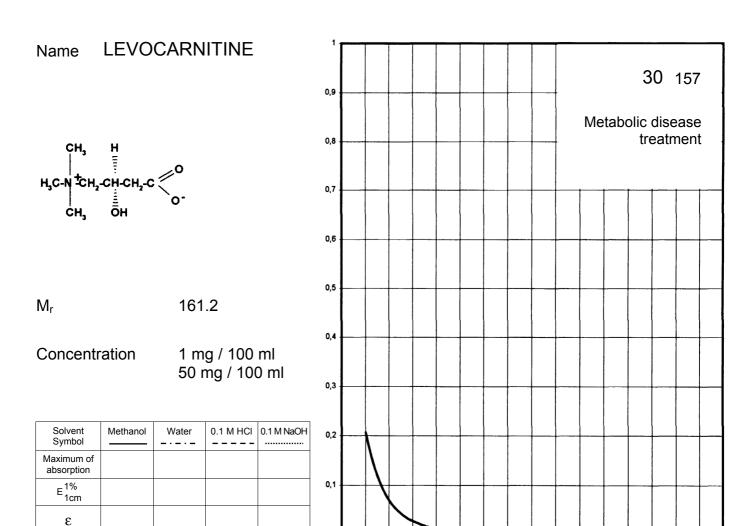

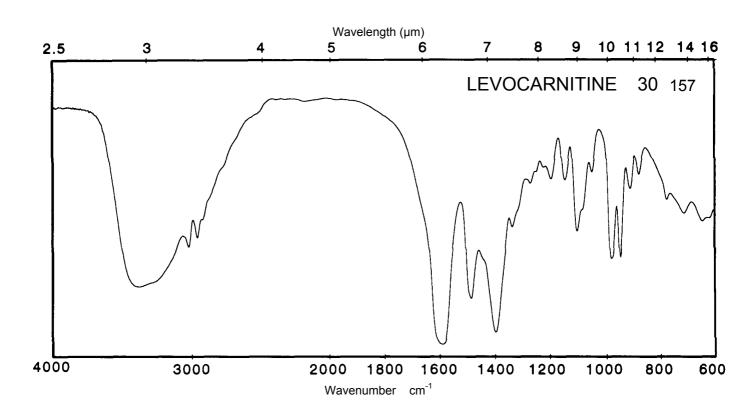



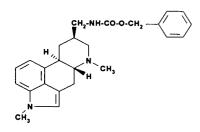



 M_{r}


403.5


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm	290 nm	289 nm	289 nm
	225 nm	225 nm	225 nm	227 nm
E ^{1%}	181	142	156	155
1cm	895	717	783	697
3	7300	5750	6300	6250
	36100	28900	31600	28100

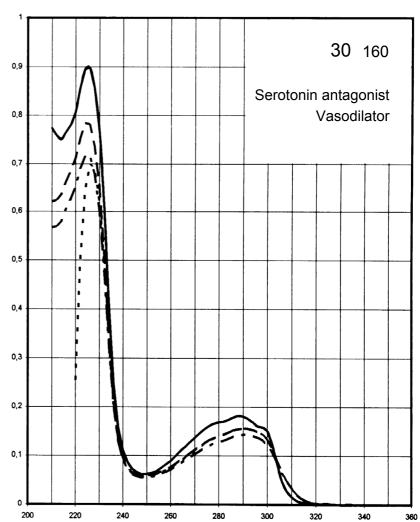


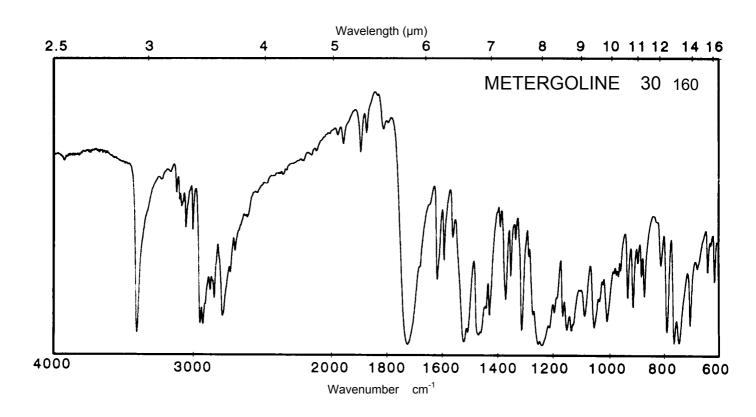

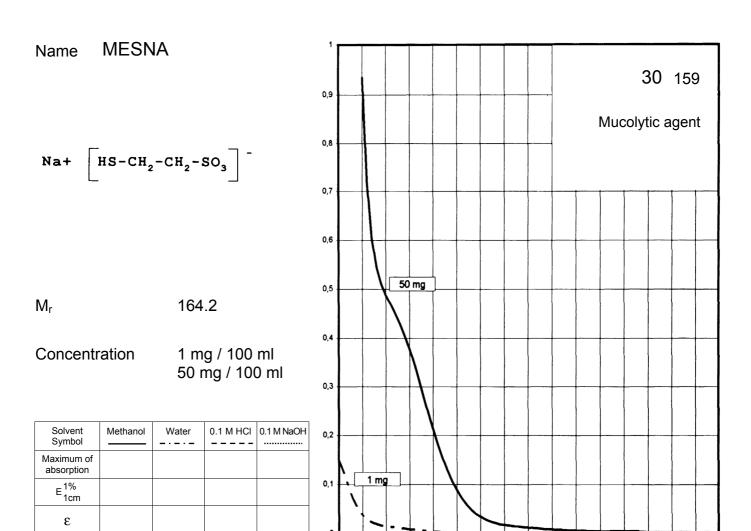

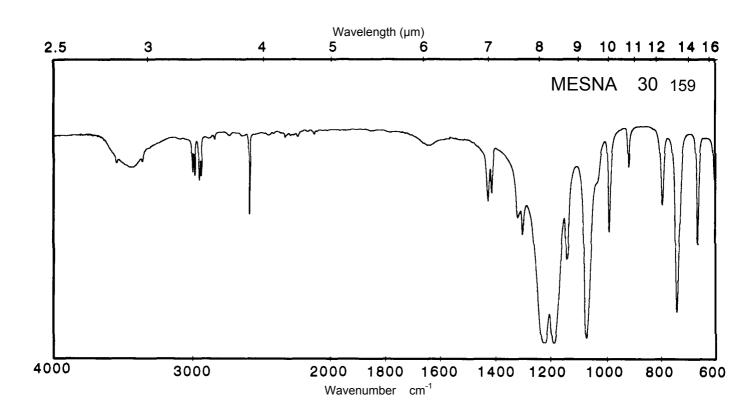




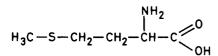


 M_{r} 149.2


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

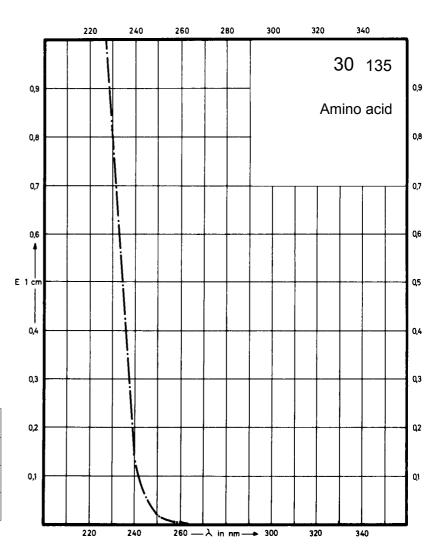

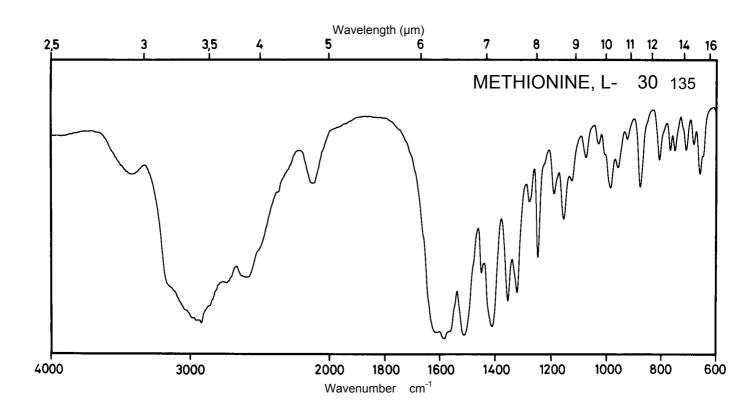



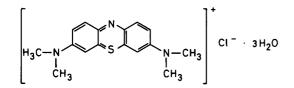

 M_{r}


373.9


Concentration

0.5 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		664 nm 292 nm		
E ^{1%} *		2400 1350		
£*		76800 43200		

^{*} Calculated on dried substance

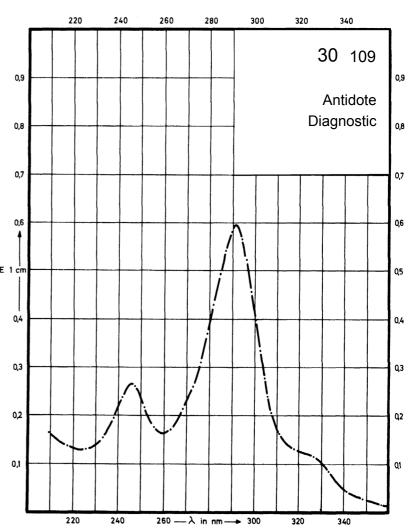

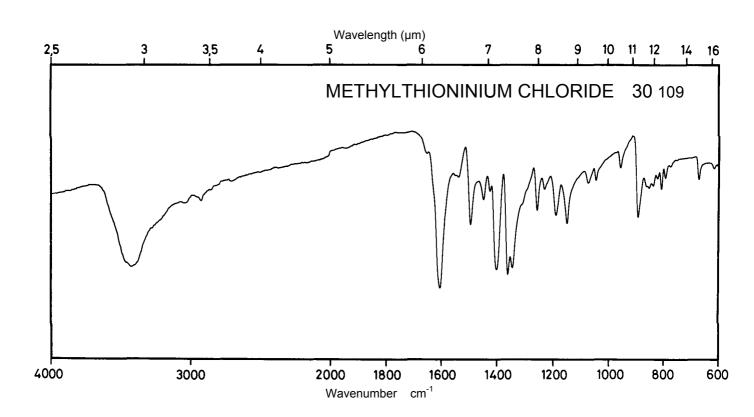

Name METHYLERGO-METRINE MALEATE

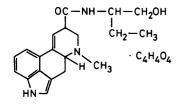



 M_{r}

455.5

Concentration 1.7 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	312 nm		312 nm	309 nm
E ^{1%} 1cm	196		182	199
3	8940		8300	9050

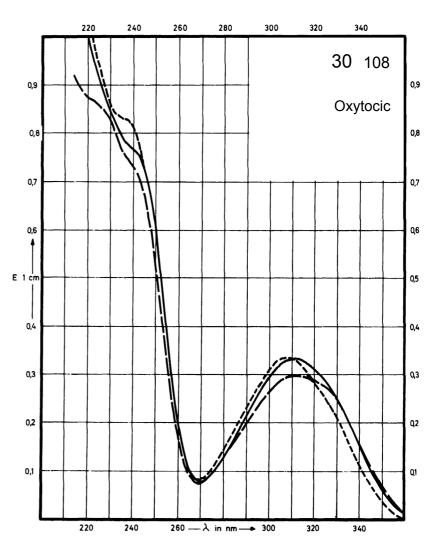

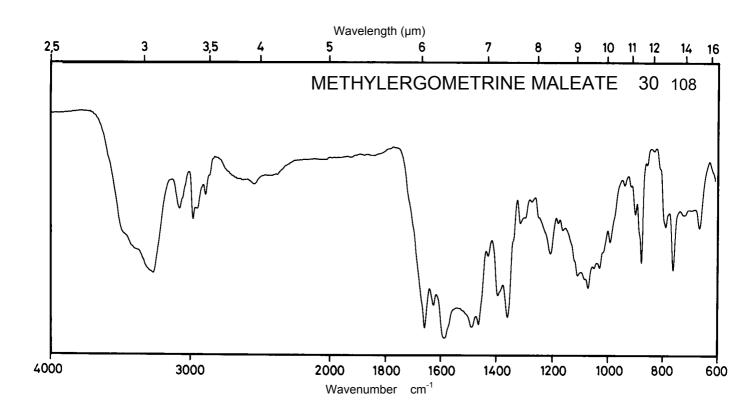



Name METAMFEPRAMONE HYDROCHLORIDE


 M_r 213.7

Concentration 1 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	249 nm		252 nm	
E ^{1%} 1cm	624		620	
ε	13330		13250	

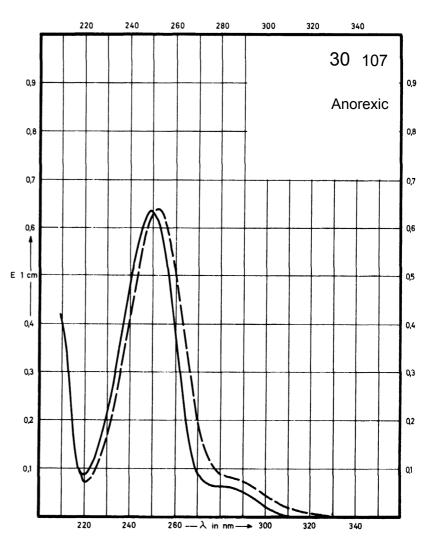


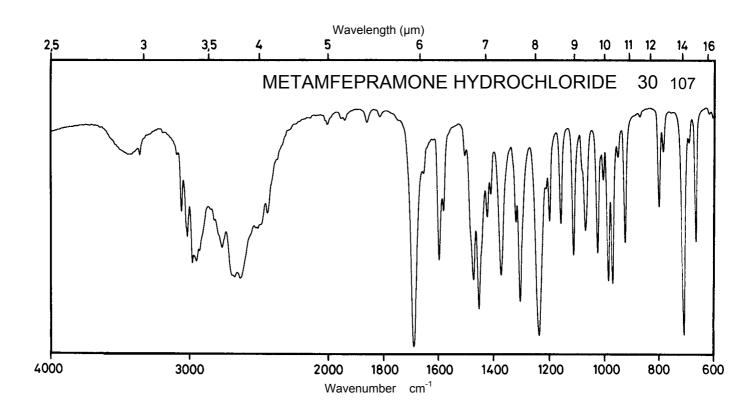

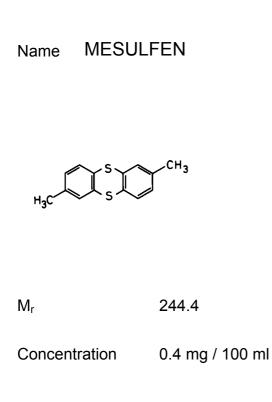


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm			
E ^{1%} 1cm	1400			
ε	34200			

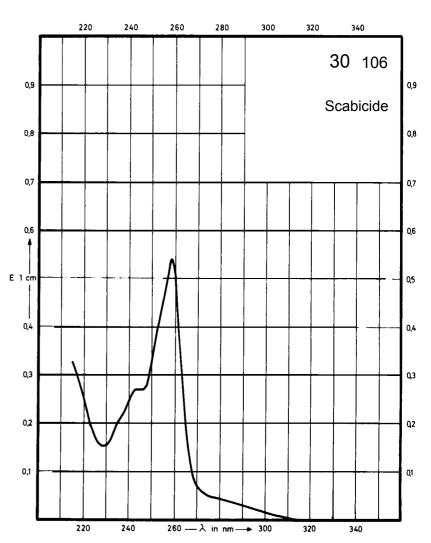

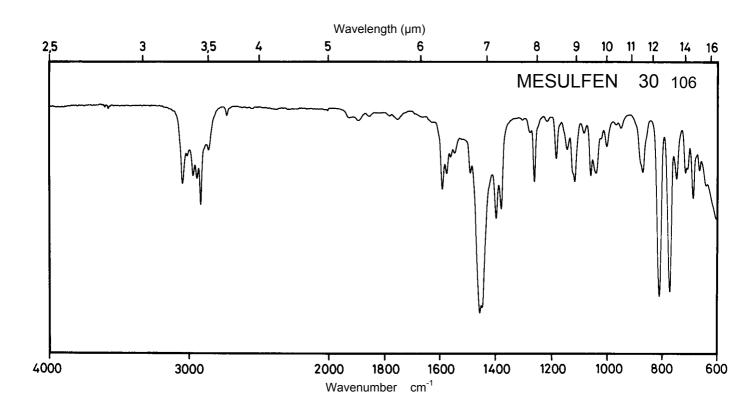

Name MEFLOQUINE HYDROCHLORIDE

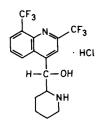



 M_{r}


414.8


Concentration 0.5 mg / 100 ml 2.5 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 282 nm 222 nm		317 nm 283 nm 222 nm	316 nm 291 nm 221 nm
E ^{1%} 1cm	67 140 1095		81 137 1010	87 143 1126
ε	2760 5880 45450		3350 5700 41820	3590 5940 46700

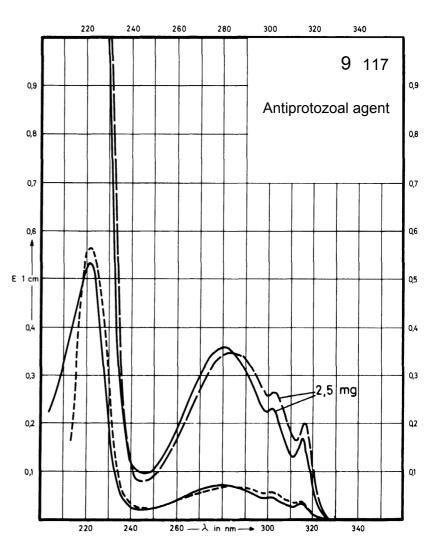


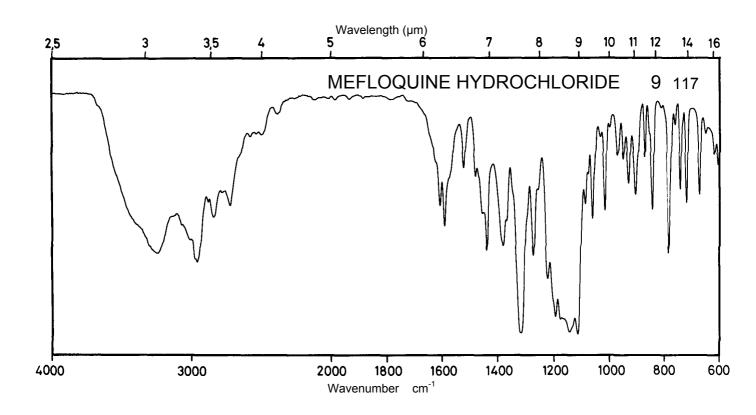

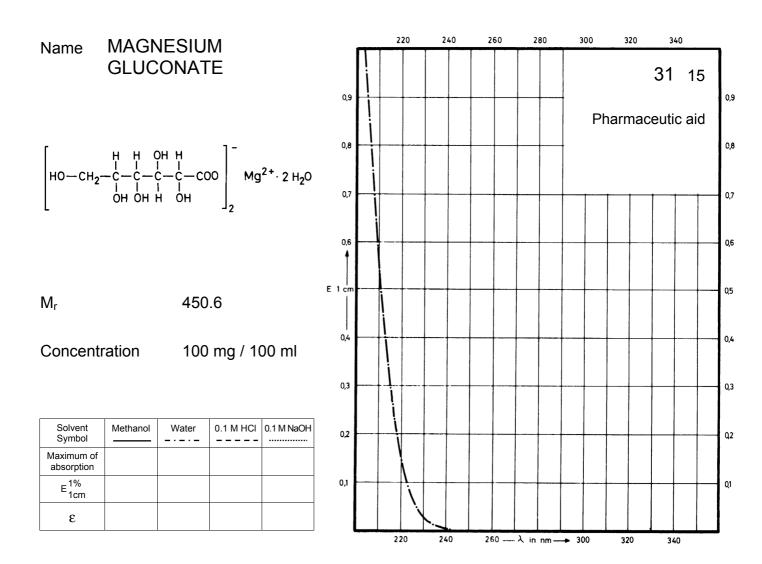


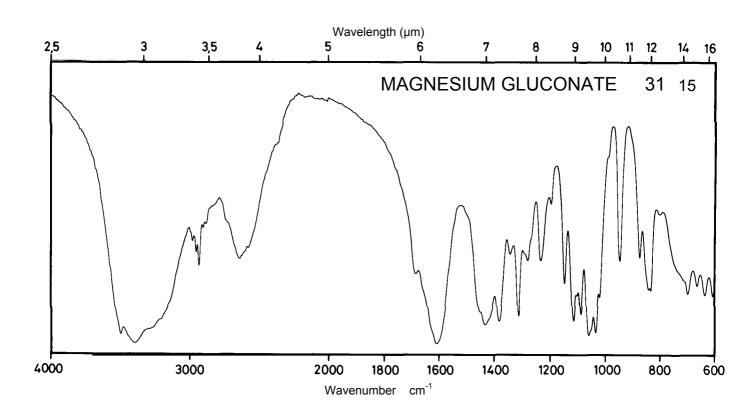

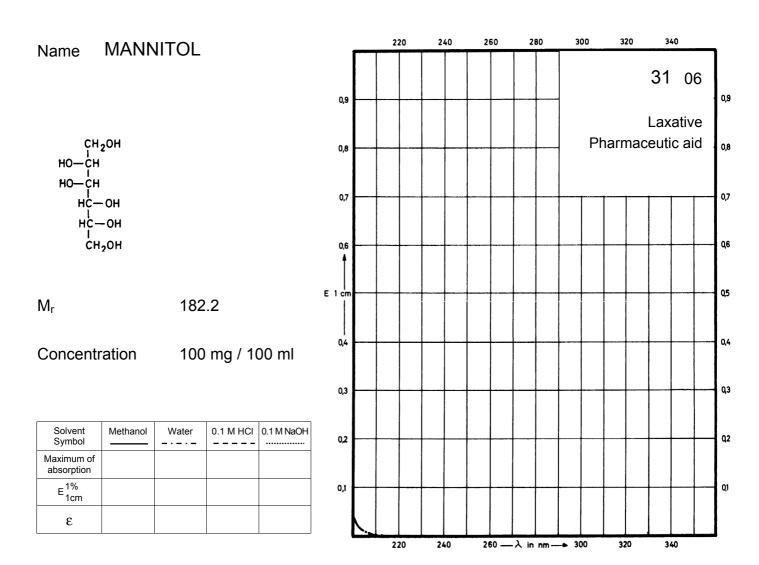

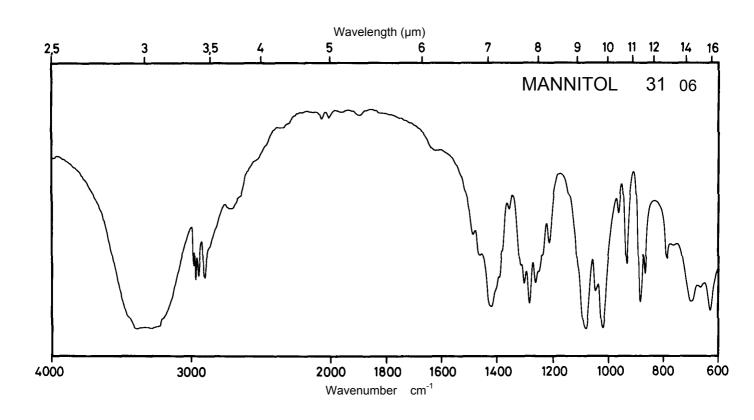



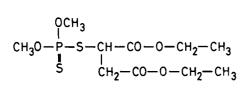



 M_{r}


330.4


Concentration 10 mg / 100 ml 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

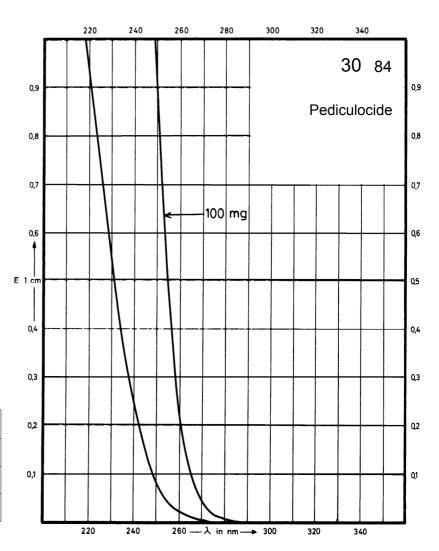

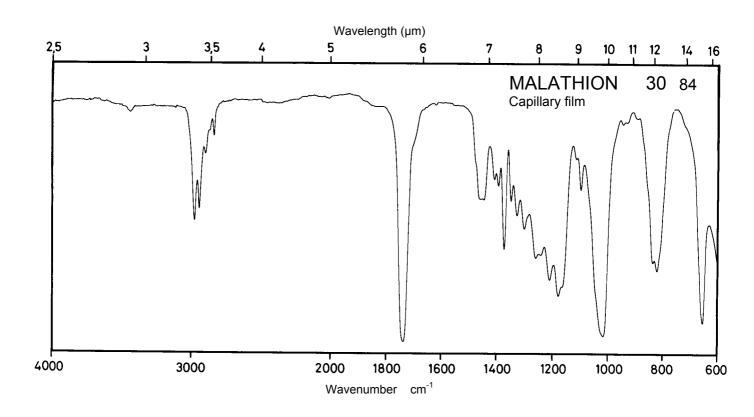

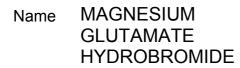


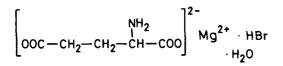


 M_{r} 268.4

Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

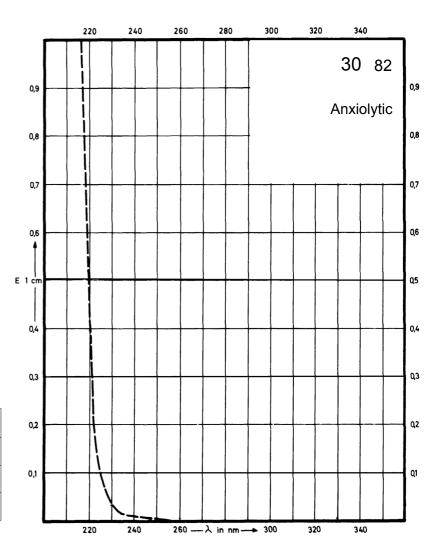

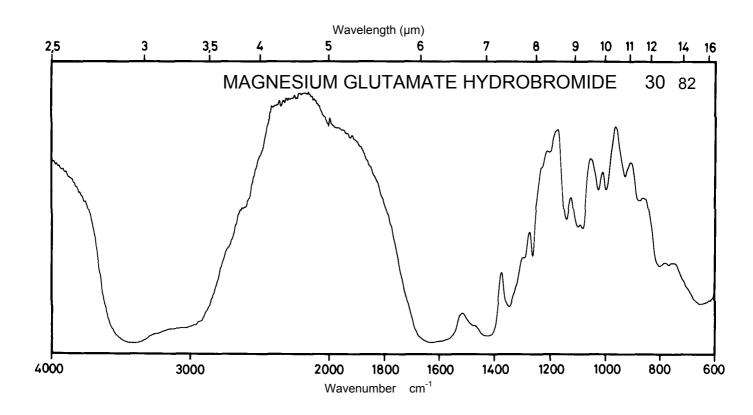


Name METHYSERGIDE MALEATE


 M_r

Concentration 1 mg / 100 ml


469.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	322 nm	323 nm	323 nm 230 nm	320 nm 243 nm
E ^{1%} 1cm	185	166	165 530	181 510
ε	8700	7800	7700 24900	8500 23900

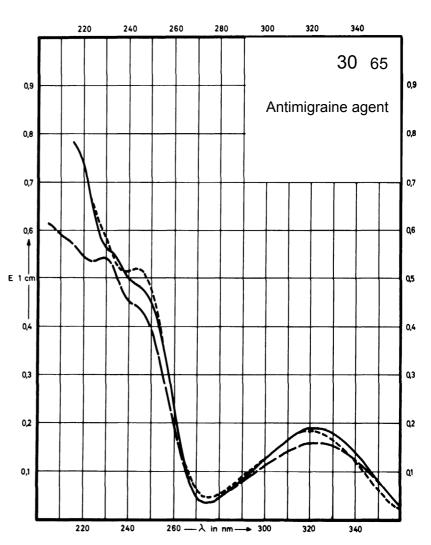

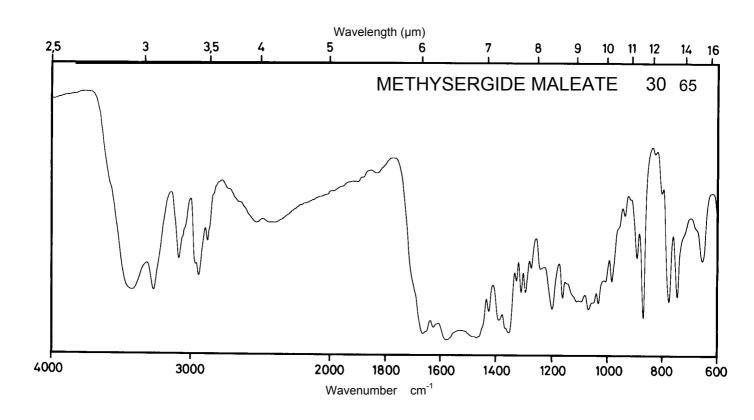


Name METIXENE HYDROCHLORIDE

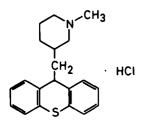


 M_r 345.9

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm	268 nm	269 nm	
E ^{1%} 1cm	326	289	291	
3	11270	9990	10060	

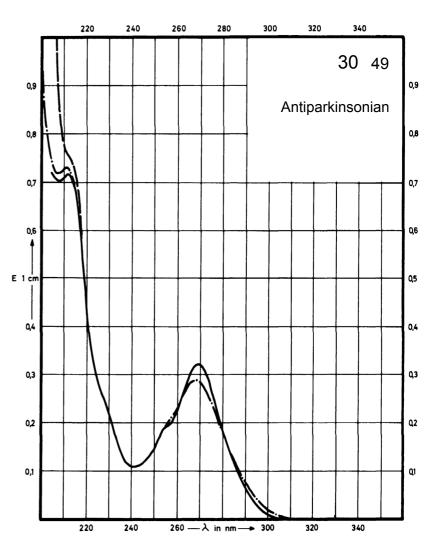

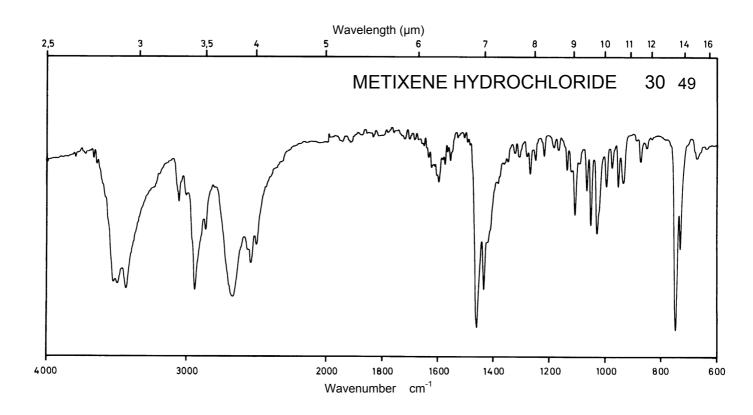


Name MEFENOREX HYDROCHLORIDE


 M_r

248.2


Concentration 75 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 253 nm	264 nm 258 nm 253 nm		
E ^{1%} 1cm	5.68 7.25 5.65	5.70 7.41 5.93		
ε	141 180	141 184 147		

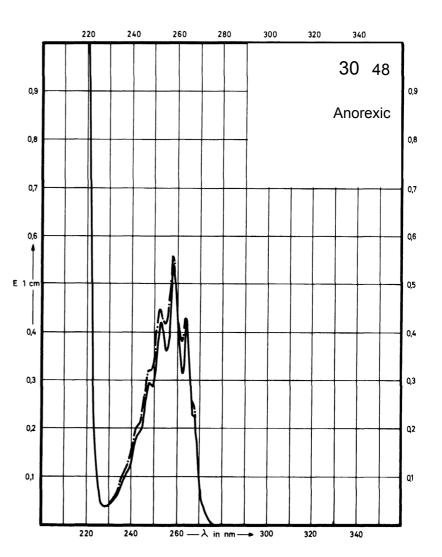

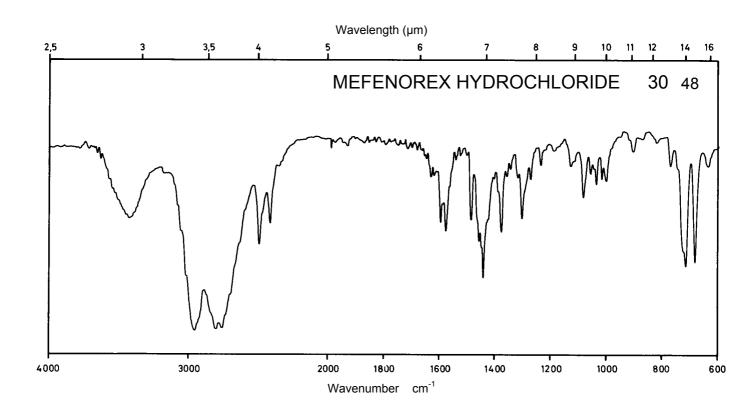




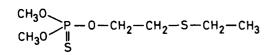


 M_r 230.3


Concentration 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

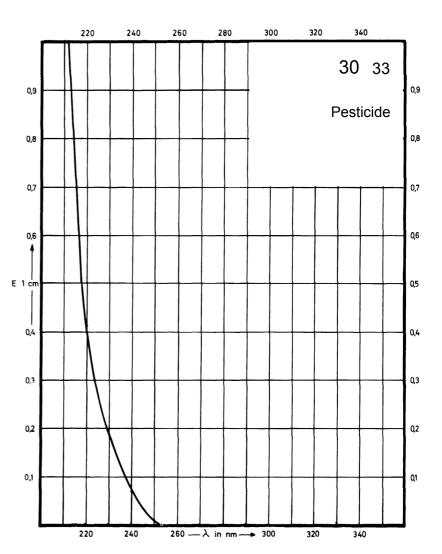

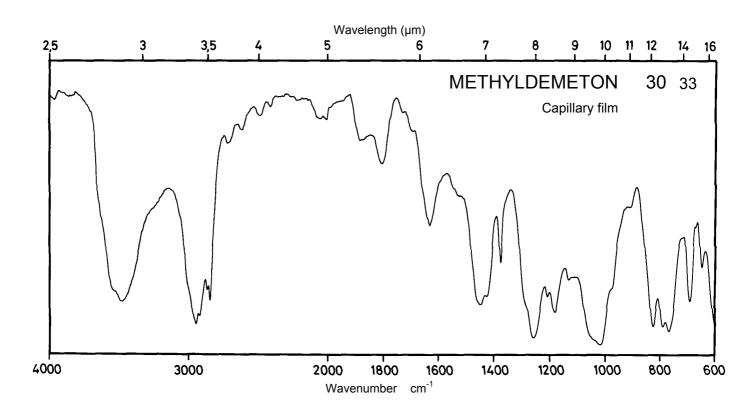



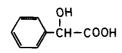



 M_{r}


152.1


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm	261 nm	262 nm	263 nm
	258 nm	257 nm	257 nm	257 nm
	252 nm	251 nm	251 nm	251 nm
E ^{1%} 1cm	9.9 12.9 10.5	10.4 13.6 10.6	11.2 13.7 10.9	9.8 13.0 10.4
ε	150	158	170	149
	196	207	208	198
	160	161	166	158

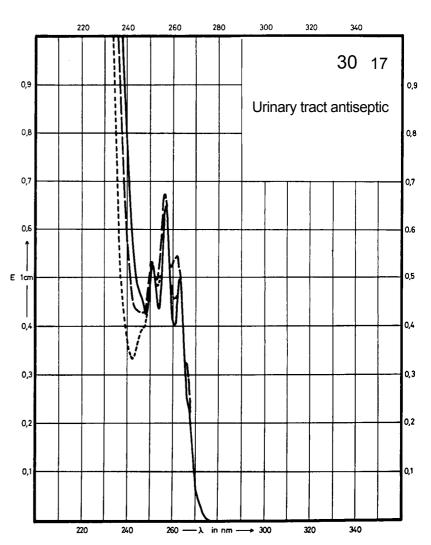


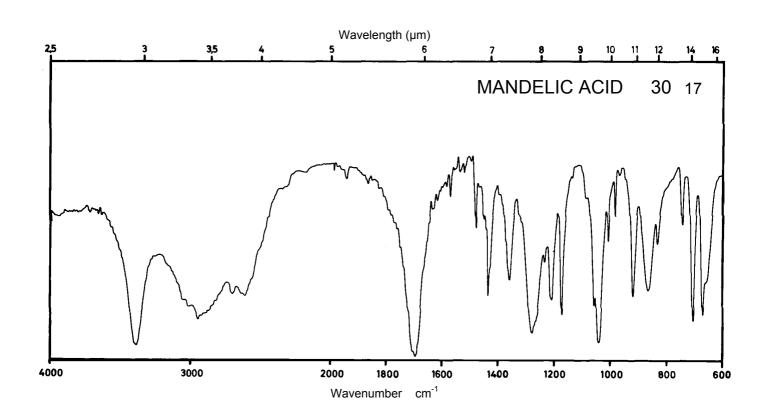

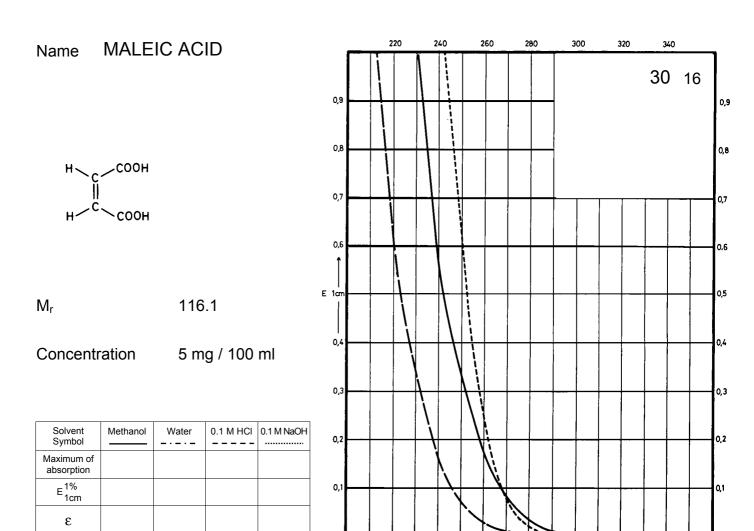

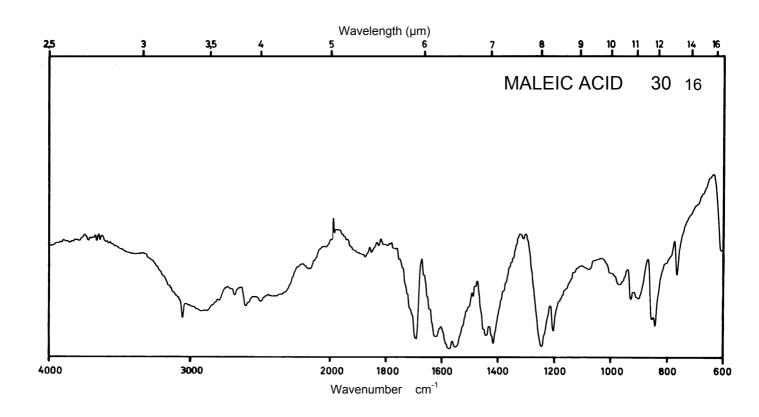



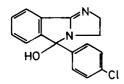



 M_r 284.7

Concentration 1 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH*
Maximum of absorption	274 nm 268 nm		271 nm	275 nm 269 nm
E ^{1%} 1cm	212 220		500	173 170
ε	6040 6260		14240	4930 4840

* 1 M NaOH + Methanol (1 + 9)

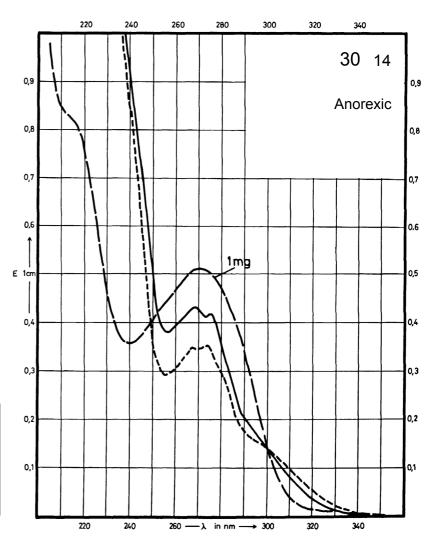

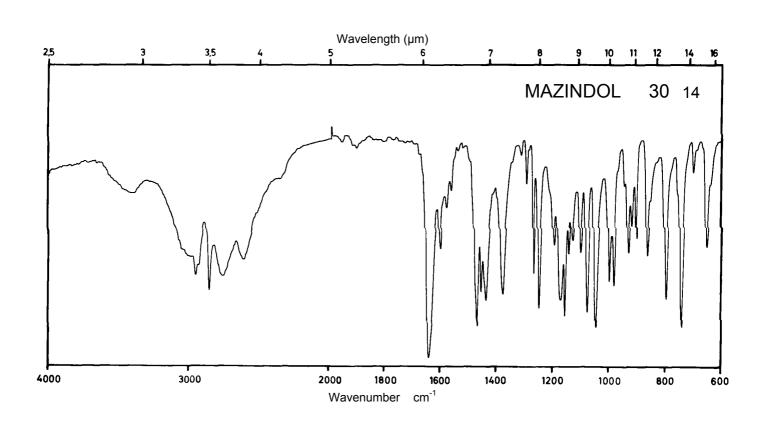


Name METHYLAMINO PHENAZONE


 M_r

217.3


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	252 nm		257 nm	268 nm
E ^{1%} 1cm	392		424	383
3	8520		9210	8320

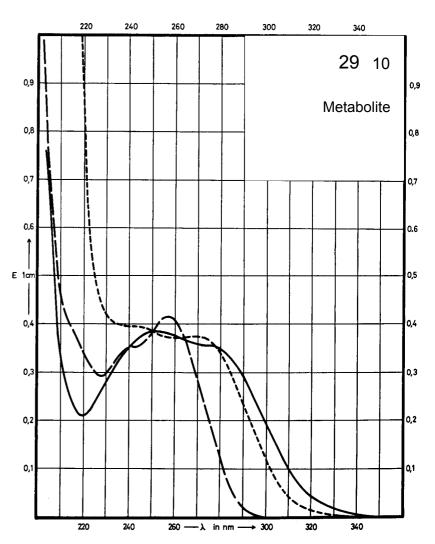

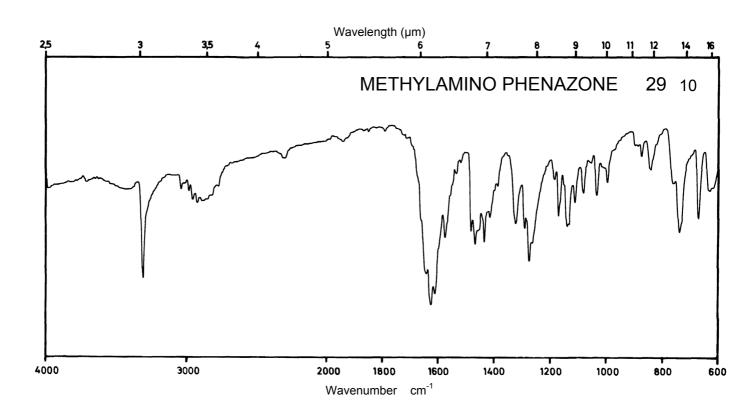



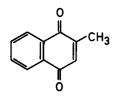


 M_{r}


172.2

Concentration


0.5 mg / 100 ml

		,	,	,
Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	330 nm 249 nm		340 nm 249 nm	
E ^{1%} 1cm	143 1118		159 1133	
ε	2460 19250		2740 19510	

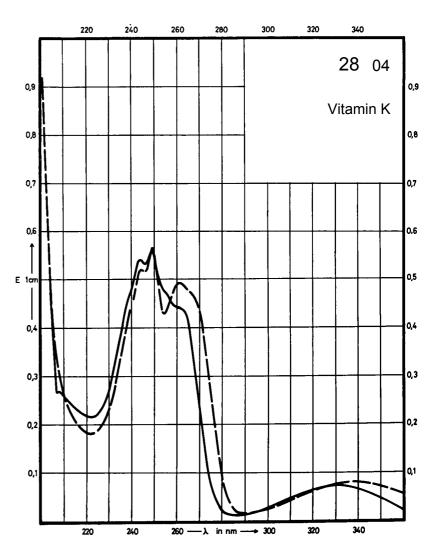

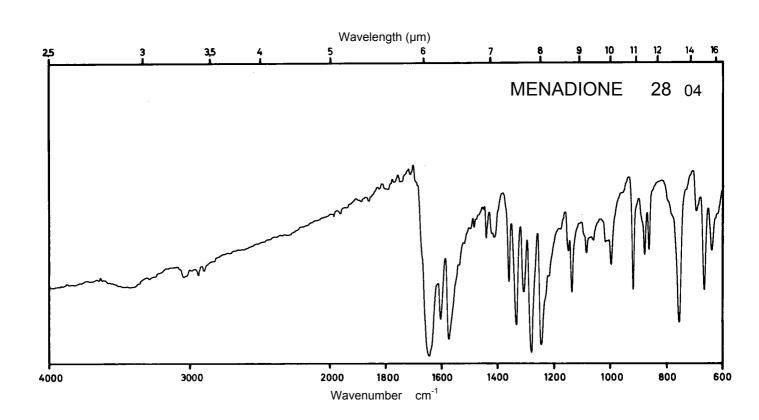




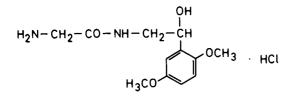


 M_{r} 290.7


Concentration 2 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 291 nm 227 nm 290 nm 290 nm absorption E^{1%} 1cm 128 264 114 114 3720 7680 3 3300 3300

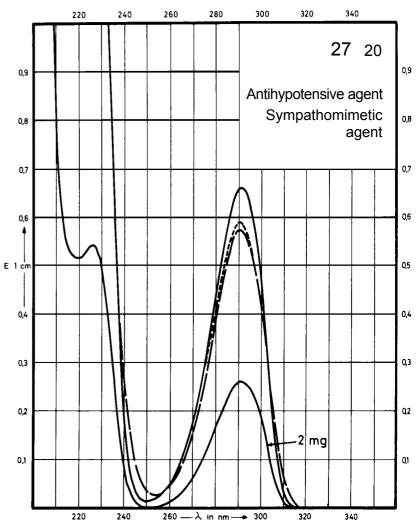

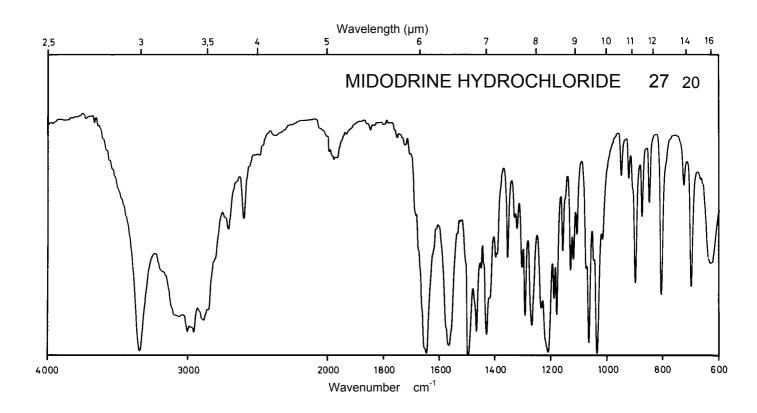

Name METARAMINOL TARTRATE

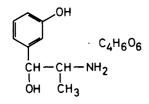



 M_r

317.3

Concentration 2.5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 218 nm		274 nm	292 nm 238 nm
E 1% 1cm	66 209		58	95 281
ε	2070 6620		1850	3000 8920

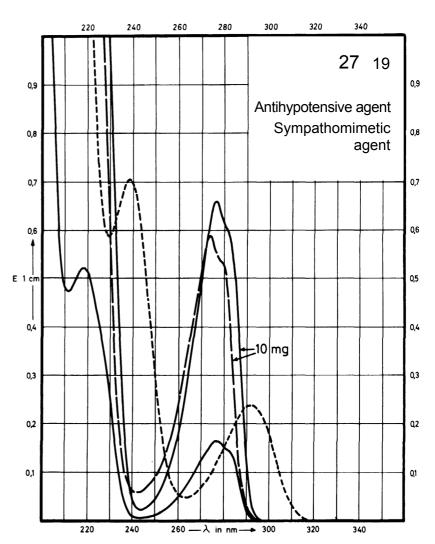

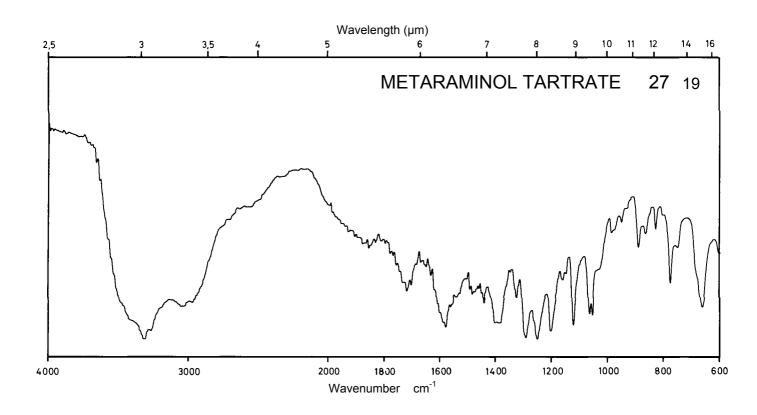


Name METHYLSYNEPHRINE HYDROCHLORIDE


 M_{r}

217.7


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 225 nm		273 nm 221 nm	290 nm 242 nm
E ^{1%} 1cm	75 465		63 414	109 634
ε	1630 10120		1370 9010	2370 13800

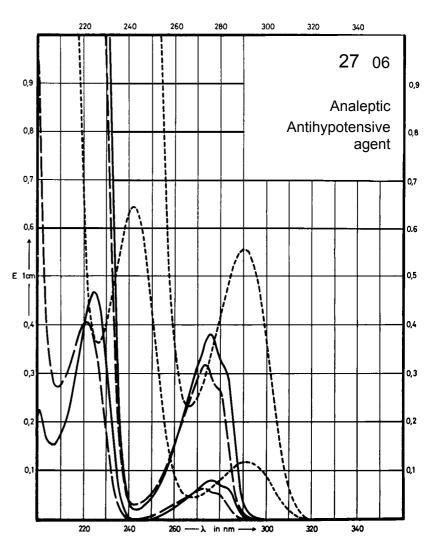

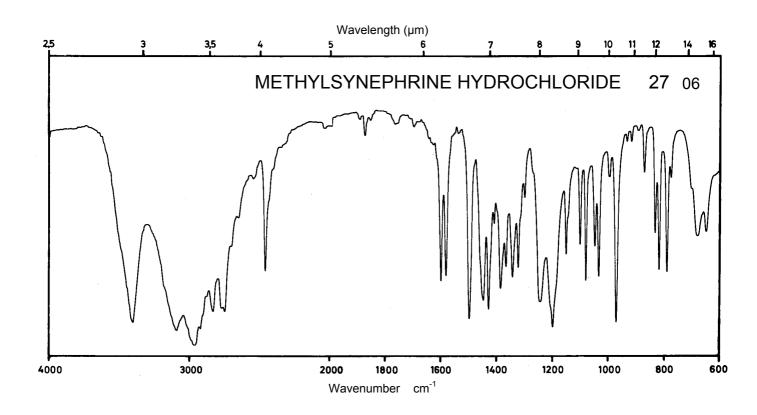

Name METACLAZEPAM HYDROCHLORIDE

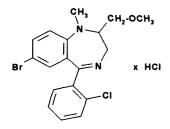



 M_r


430.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	462 nm 251 nm	460 nm 253 nm	Decom- position observed
E ^{1%} 1cm		113 513	147 574	
ε		4870 22100	6300 24700	

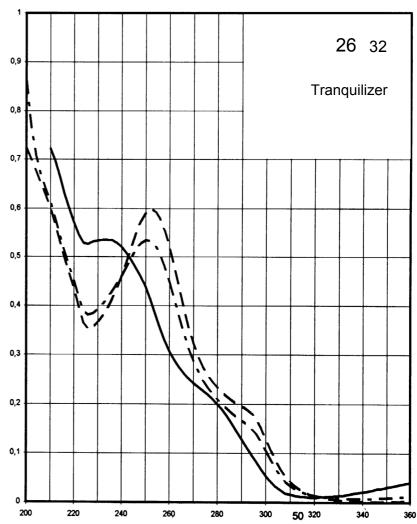

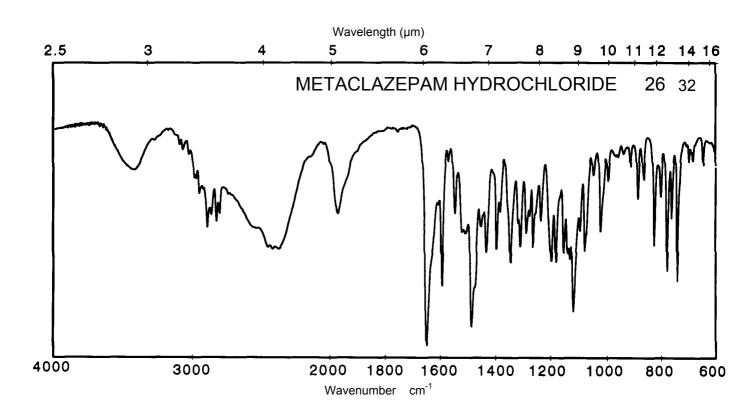


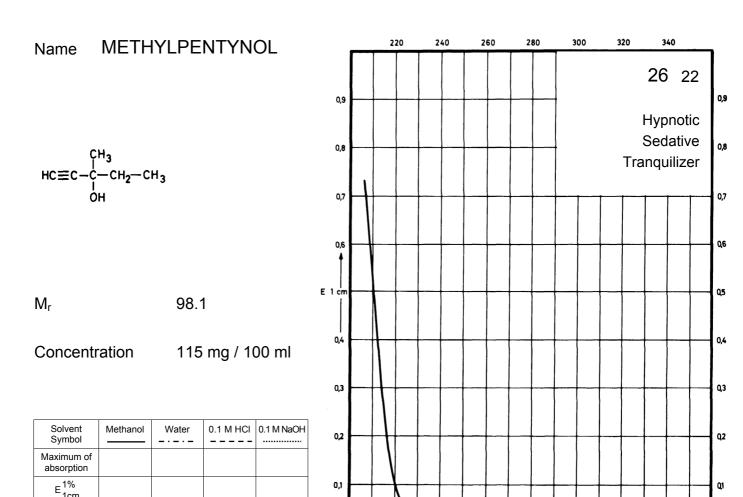

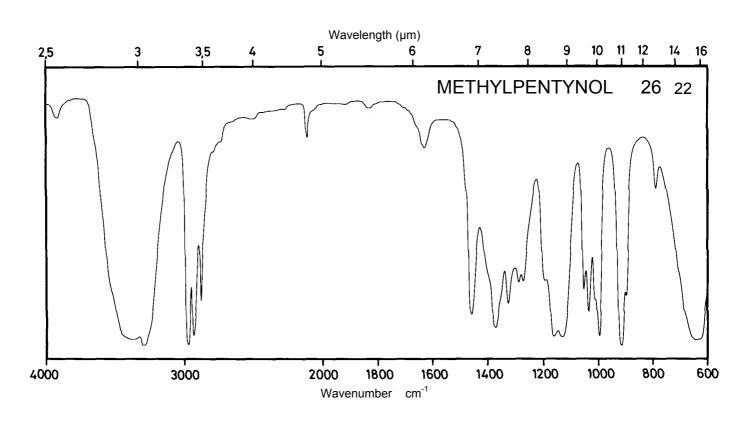

ε



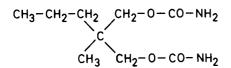



 M_r

218.3


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

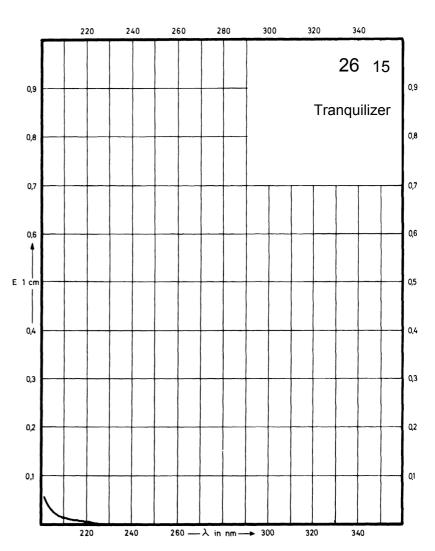

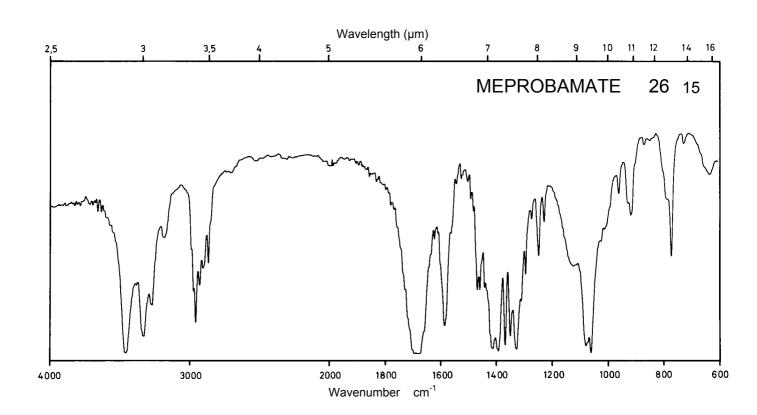




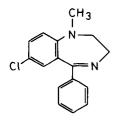


 M_r 270.8


Concentration 1 mg / 100 ml 10 mg / 100 ml

_				
Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	357 nm 231 nm		450 nm 253 nm	330 nm
E ^{1%} 1cm	69 778		163 864	56
3	1870 21060		4410 23400	1500

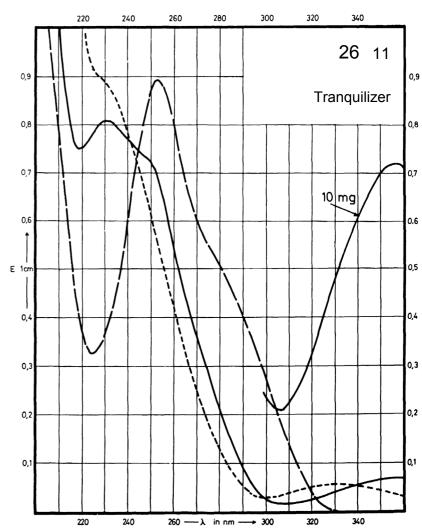

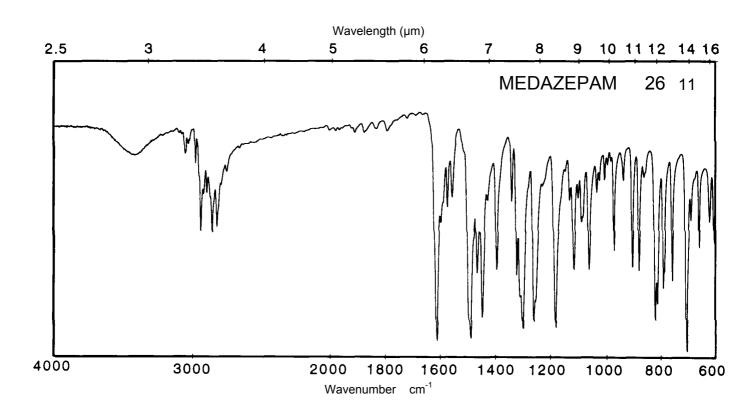



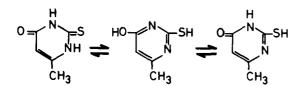



 M_r

Concentration 0.5 mg / 100 ml


142.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm		276 nm	256 nm
E ^{1%} 1cm	1150		1190	850
ε	16350		16920	12090

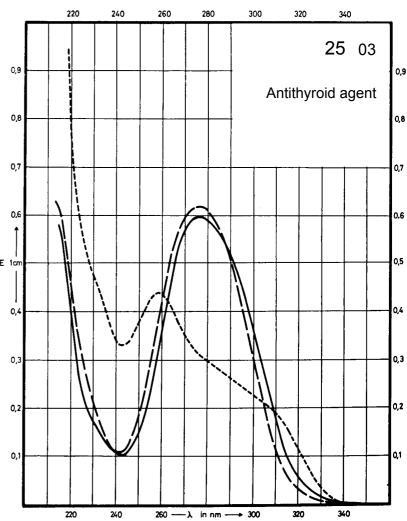

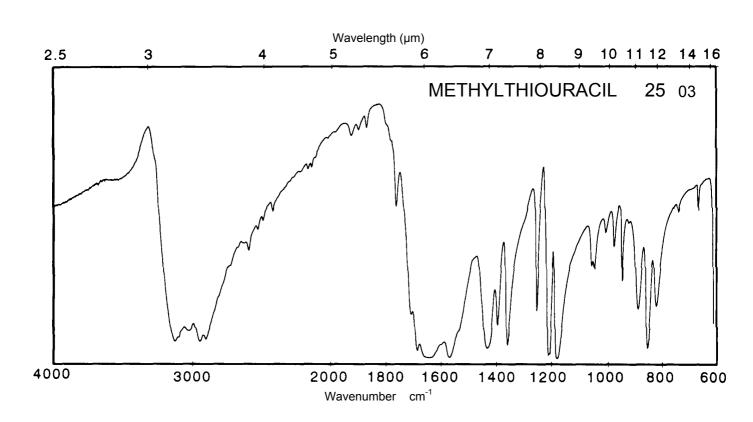




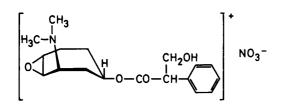



 M_r 380.4


Concentration 105 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 257 nm 251 nm	263 nm 257 nm 251 nm		
E ^{1%} 1cm	3.66 4.91 4.14	3.74 4.97 4.23		
ε	139 187	142 189		

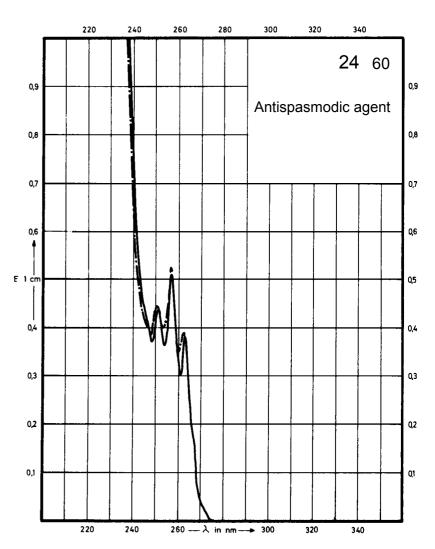

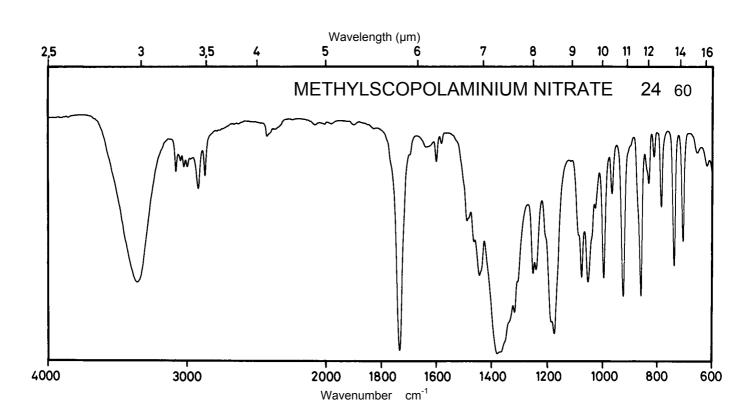

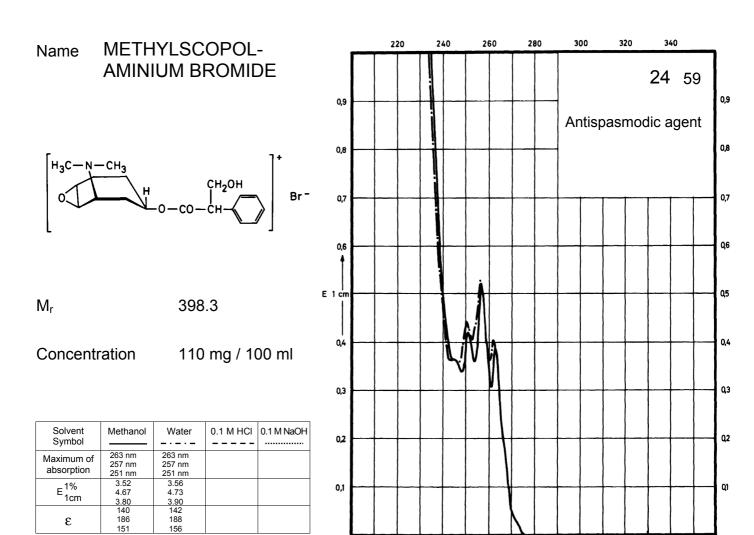

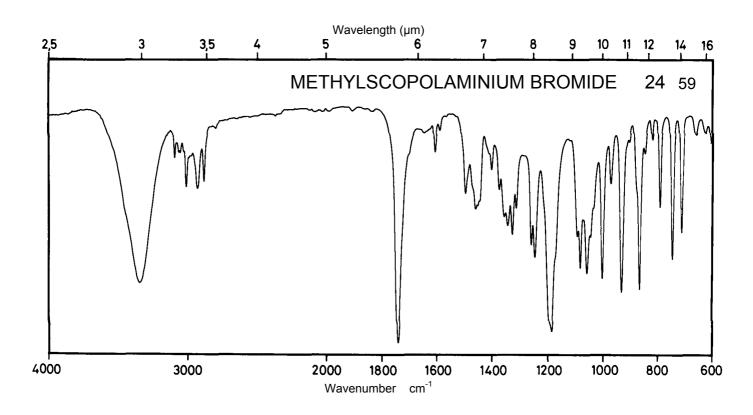

220

240

— λ in nm — 300

186 151

ε


Name MEBEVERINE HYDROCHLORIDE


 M_{r}

466.0


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm 221 nm		263 nm 220 nm	
E ^{1%} 1cm	286 695		268 621	
3	13300 32400		12500 28900	

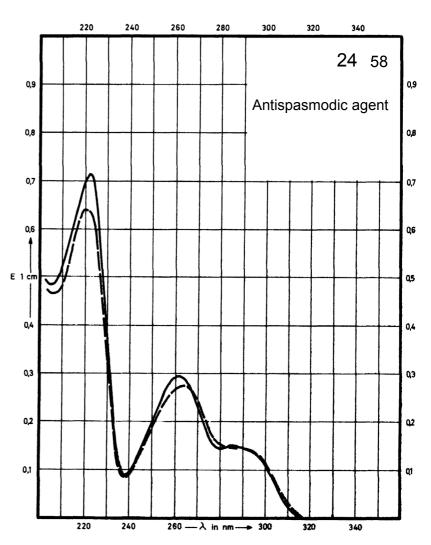

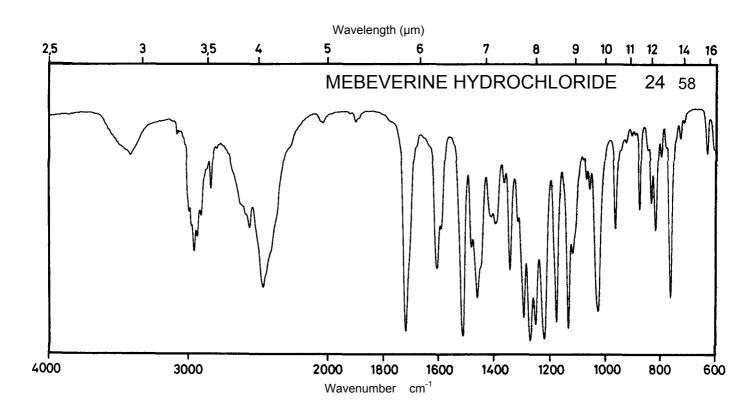




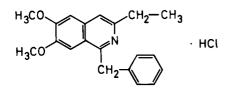


 M_r 343.9


Concentration 0.26 mg / 100 ml 2.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			252 nm	239 nm
E ^{1%} 1cm			1906	1830
ε			65500	62900

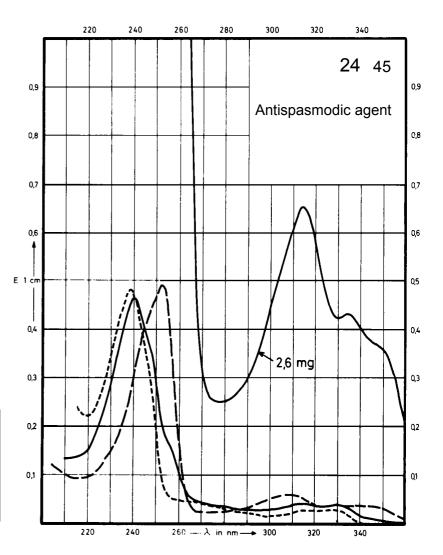

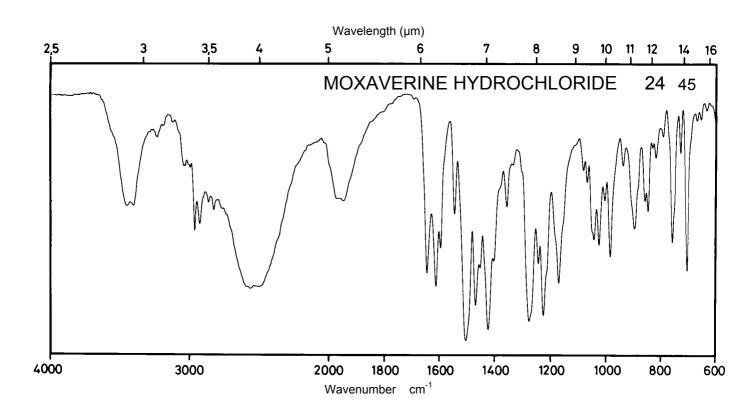




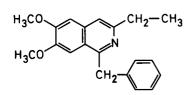


 M_r 307.4


Concentration 0.2 mg / 100 ml 2.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	332 nm 319 nm 240 nm		252 nm	239 nm
E ^{1%} 1cm	142 127 2230		2170	2100
ε	4360 3900 68600		66700	64500

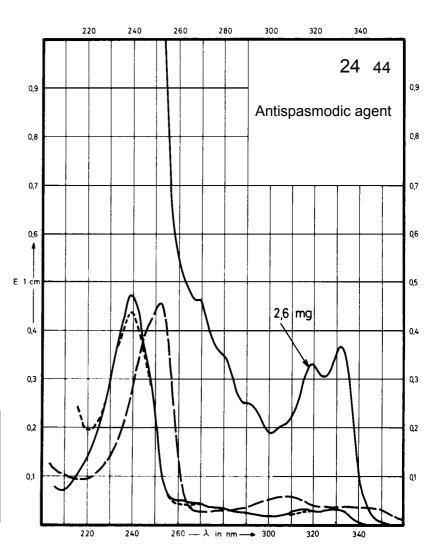

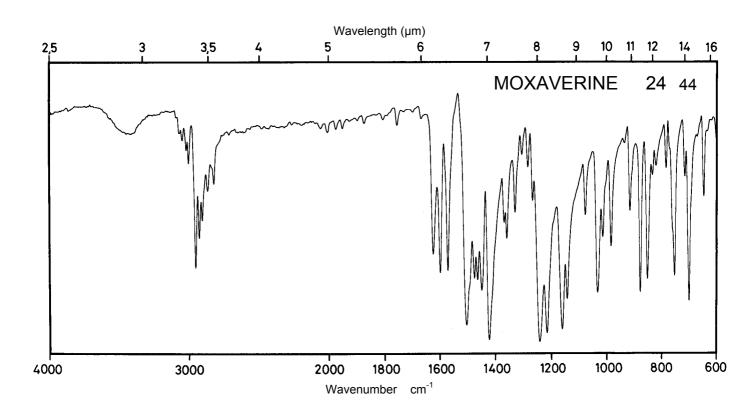




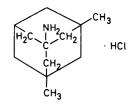


 M_r 215.8


Concentration 120 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
ε				

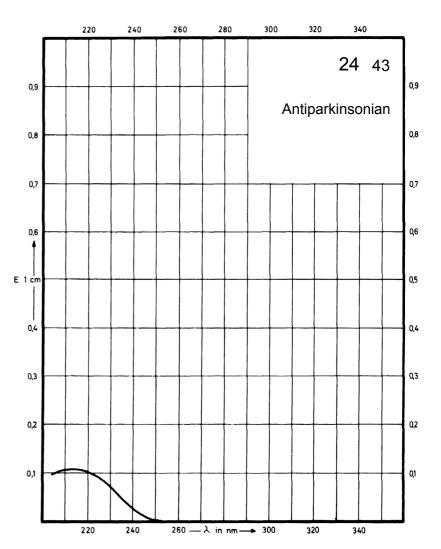

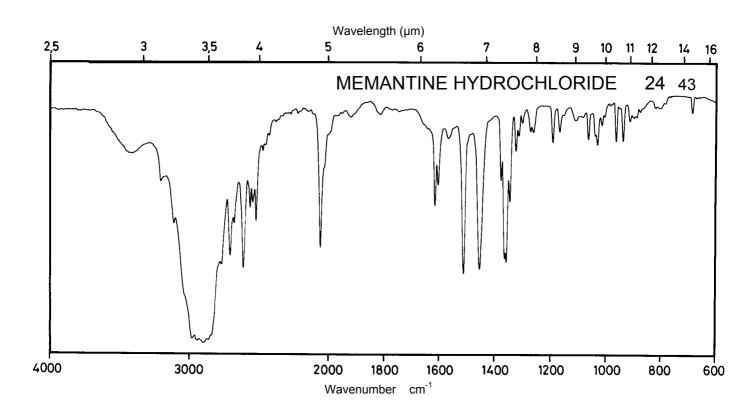



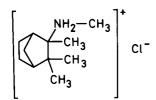



 M_r

203.8


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

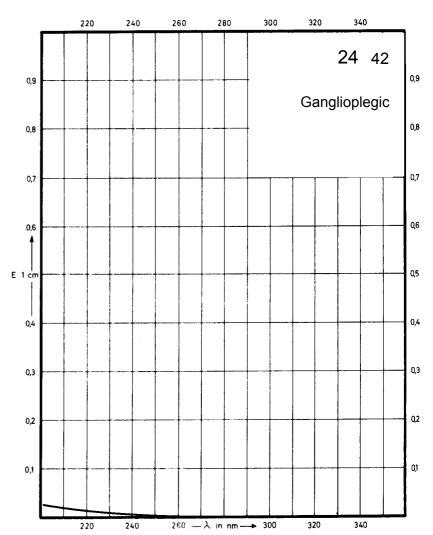

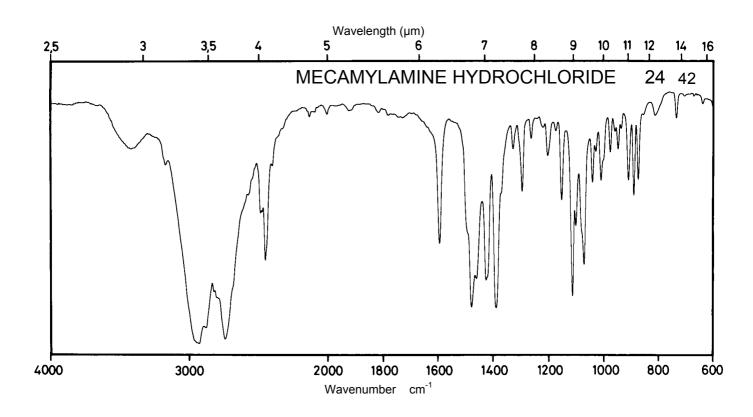



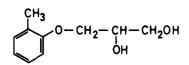



 M_{r}

182.2

Concentration 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm	276 nm	276 nm	276 nm
	272 nm	270 nm	270 nm	270 nm
E ^{1%}	86	71	71	71
1cm	95	81	81	81
ε	1570	1290	1290	1290
	1730	1480	1480	1480

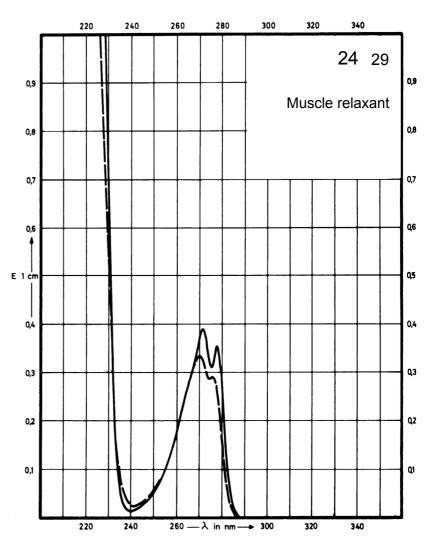

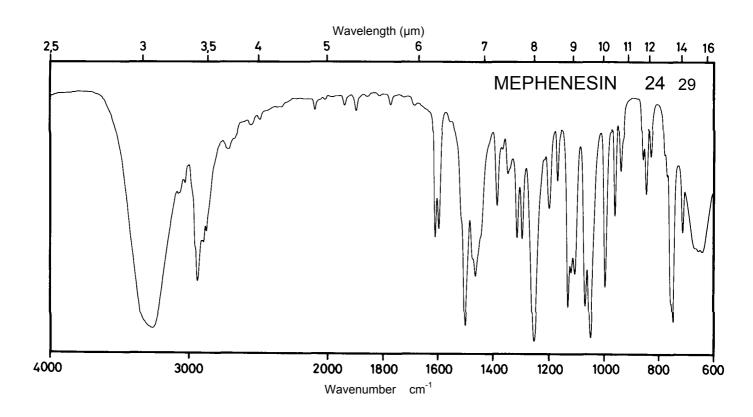



Name METHOCARBAMOL

 M_r 241.2


Concentration 2 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 225 nm		274 nm	274 nm
E ^{1%} 1cm	96 303		95	95
ε	2310 7310		2290	2290

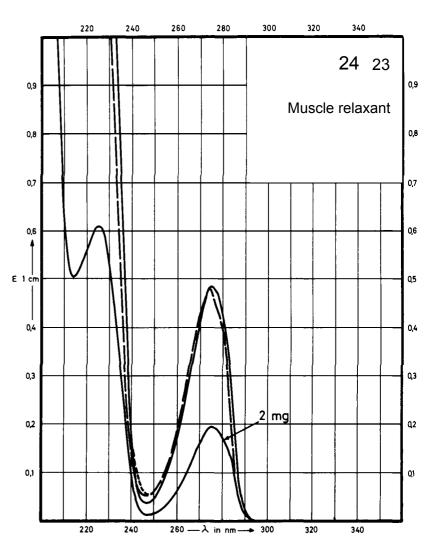

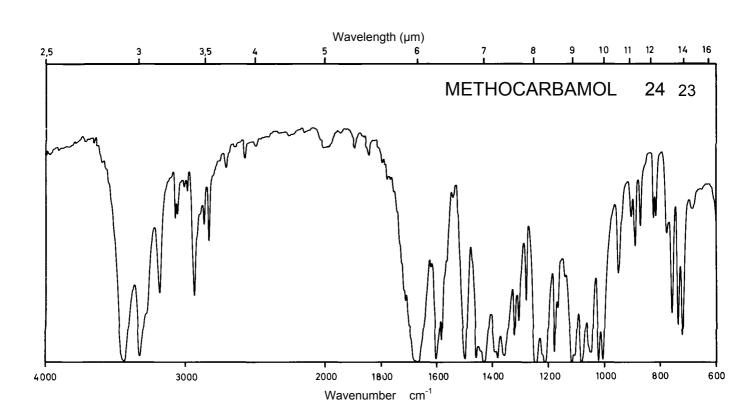




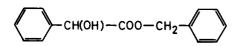


242.3


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm			
E ^{1%} 1cm	16.2			
ε	390			

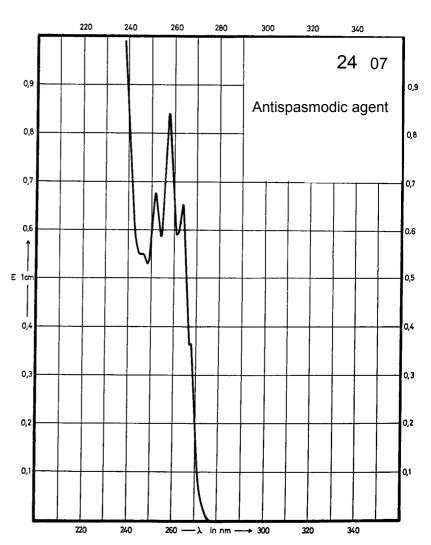

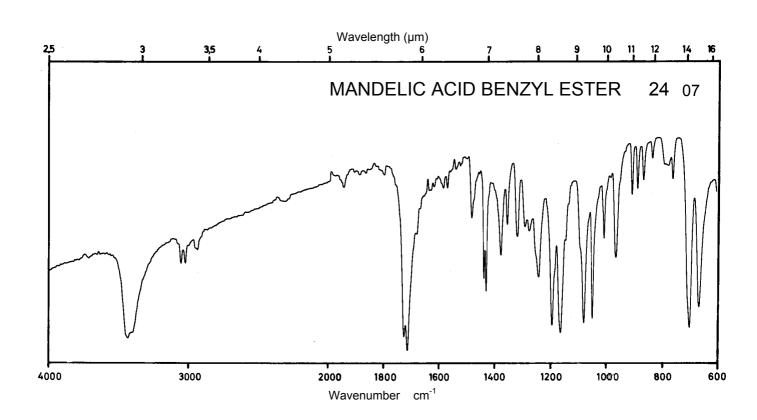




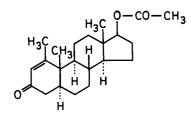


 M_r 344.5


Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm			
E ^{1%} 1cm	387			
3	13300			

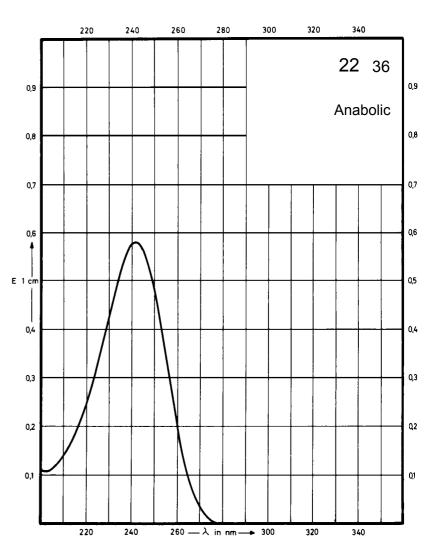

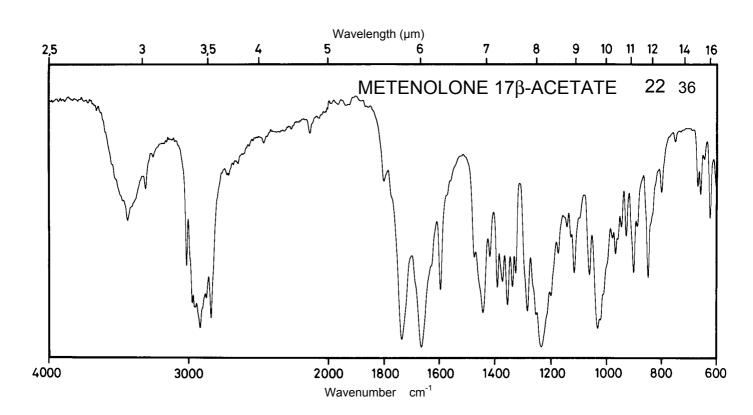




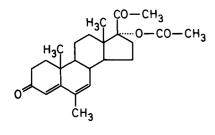


384.5


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm	298 nm		
E ^{1%} 1cm	631	605		
ε	24300	23300		

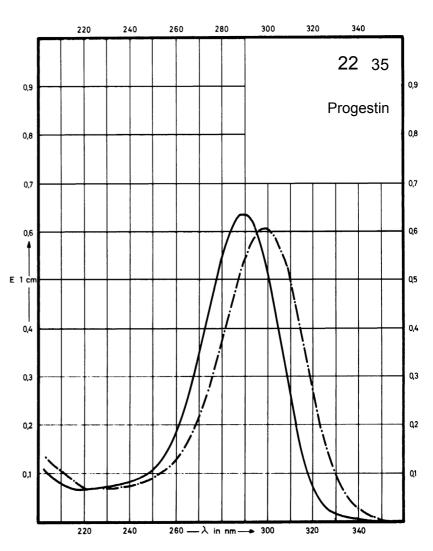

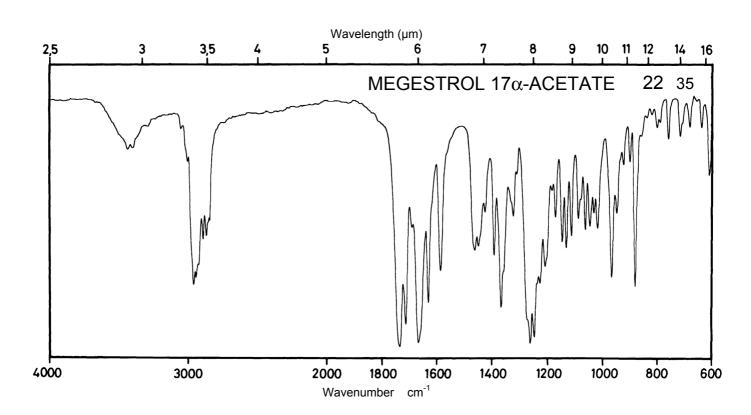




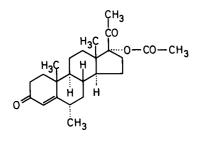


386.5


Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm			
E ^{1%} 1cm	416			
3	16100			

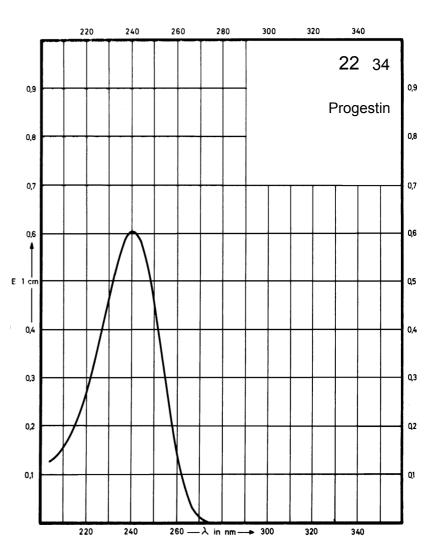

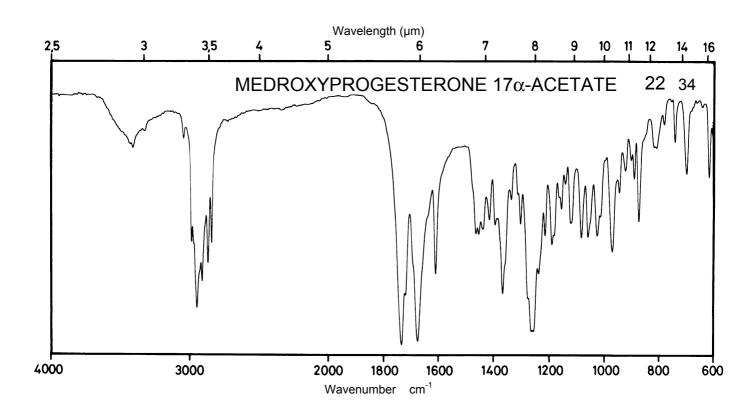




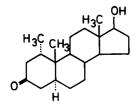


 M_r 304.5


Concentration 500 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm			
E ^{1%} 1cm	0.68			
ε	20.6			

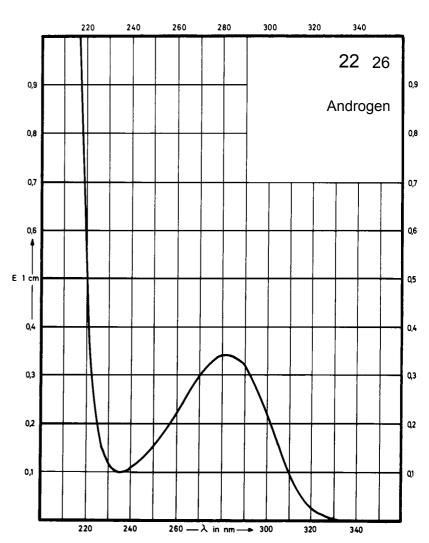

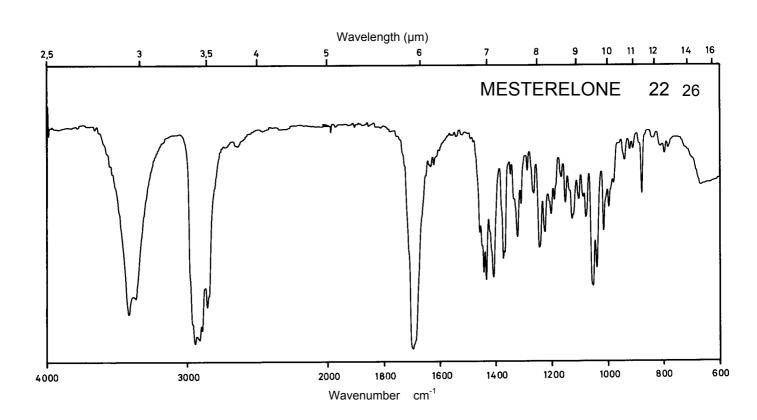




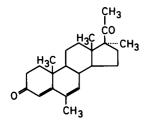


 M_r 340.5


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm			
E ^{1%} 1cm	700			
3	23830			

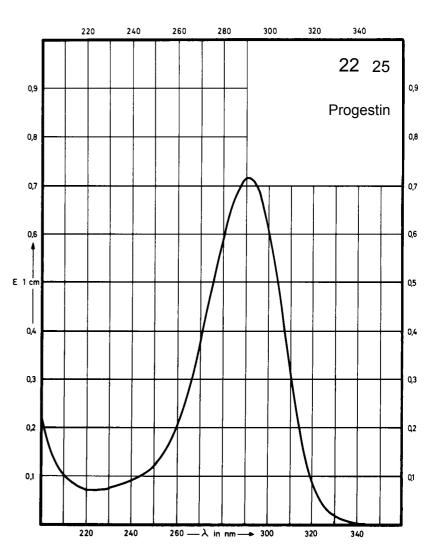

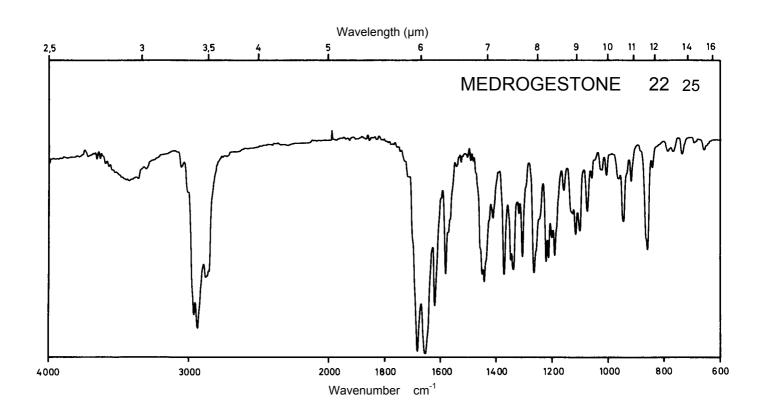




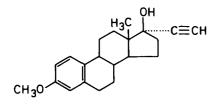


310.4


Concentration 10 mg / 100 ml

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of absorption	287 nm 278 nm			
E ^{1%} 1cm	61 65			
ε	1900 2020			

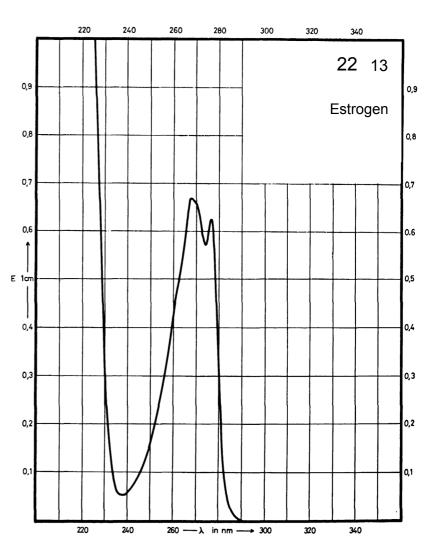

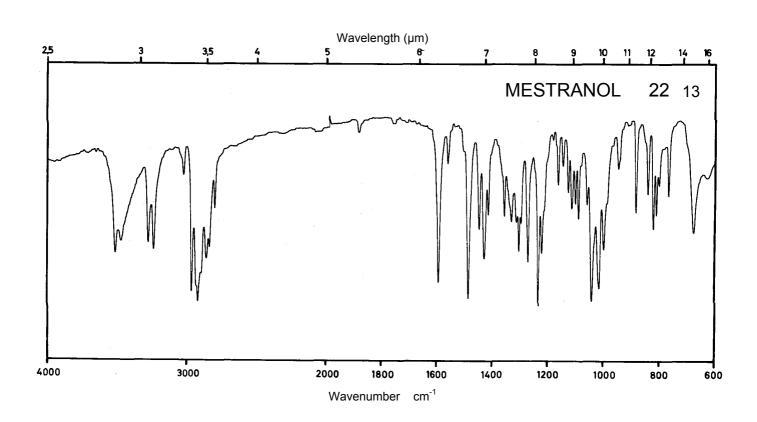




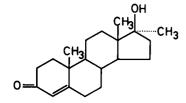


302.4


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm			
E ^{1%} 1cm	533			
ε	16120			

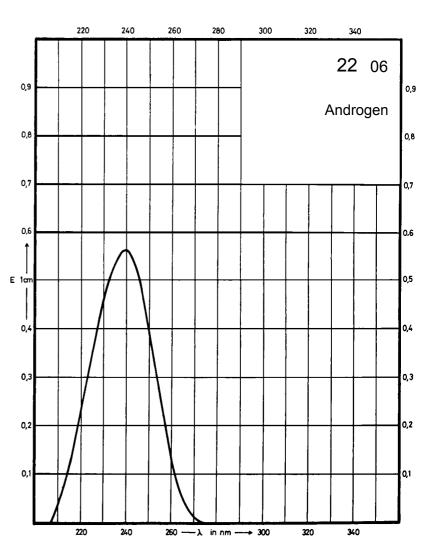

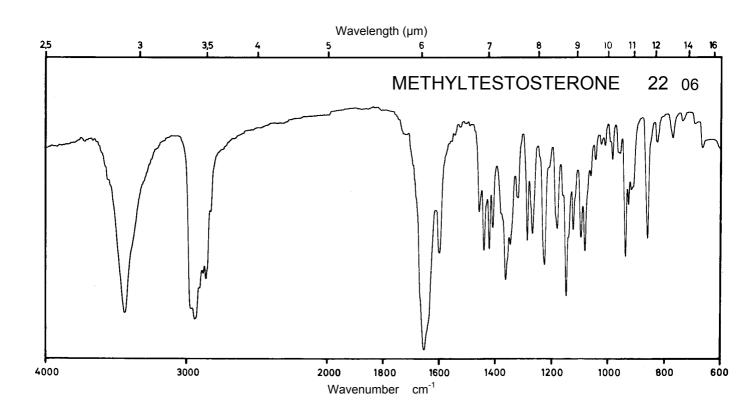



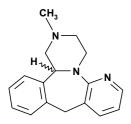



 M_r

265.4


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	293 nm	293 nm	316 nm 253 nm	294 nm
E ^{1%} 1cm	213	208	362 154	208
3	5660	5530	9610 4100	5520

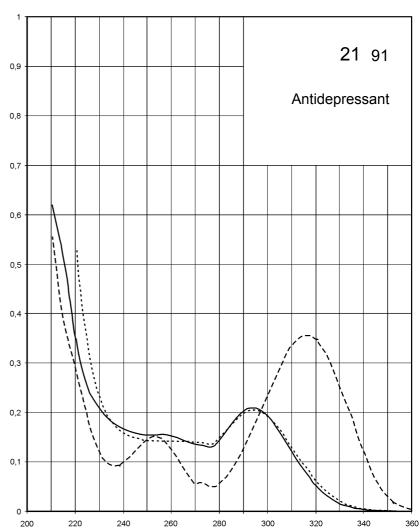

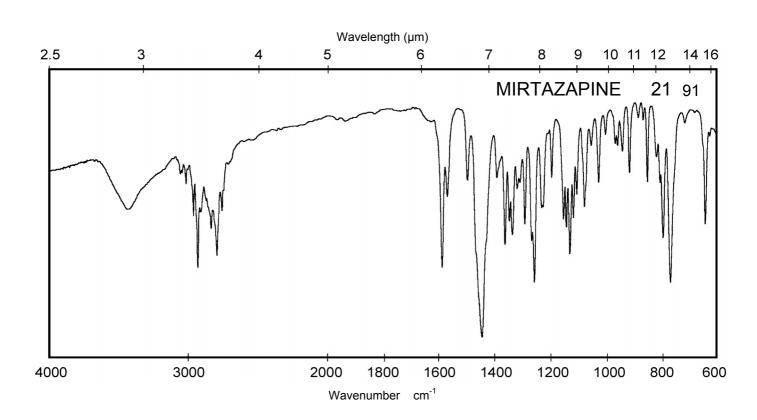




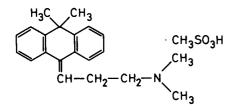


 M_r 387.5


Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm		257 nm	
E ^{1%} _{1cm}	367		348	
ε	14200		13500	

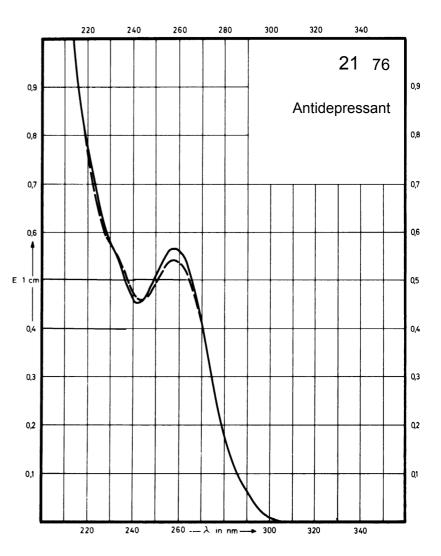

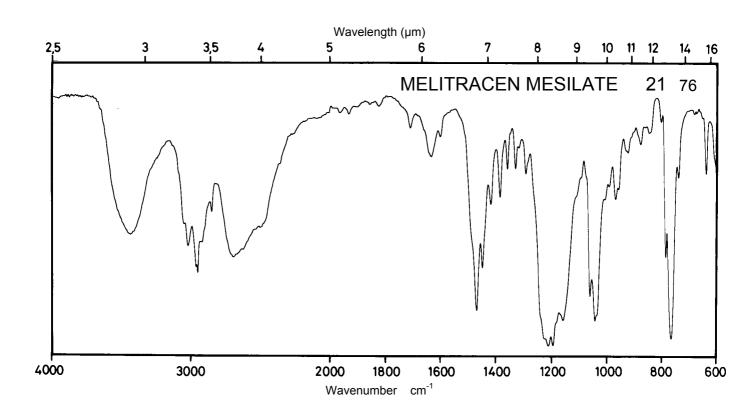


Name MELITRACEN HYDROCHLORIDE

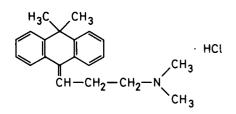



 M_r 327.9


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm		257 nm	
E ^{1%} 1cm	446		421	
ε	14600		13800	

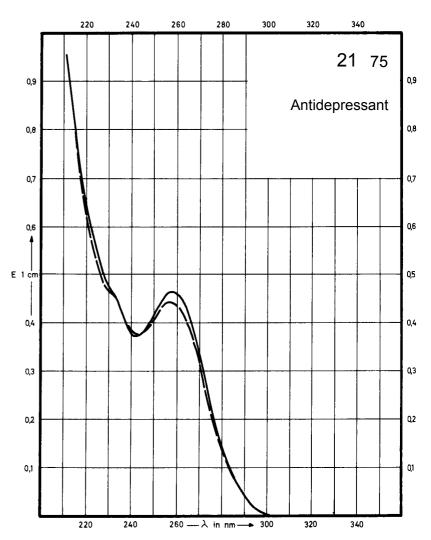

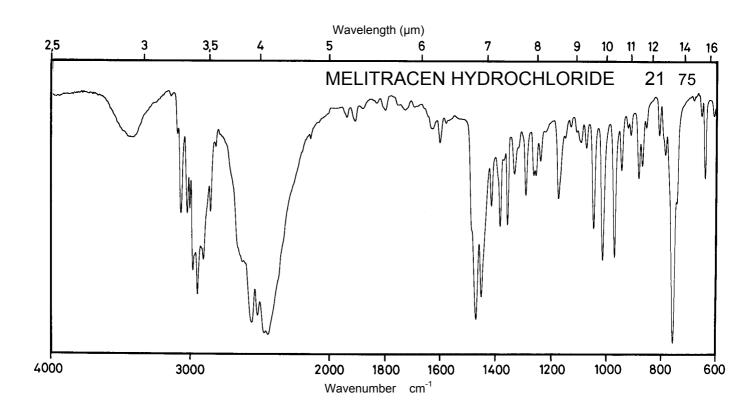

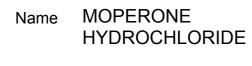


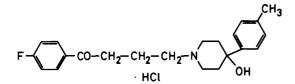


 M_{r} 391.9

Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	244 nm	248 nm	248 nm	
E ^{1%} 1cm	335	330	330	
ε	13100	12950	12950	

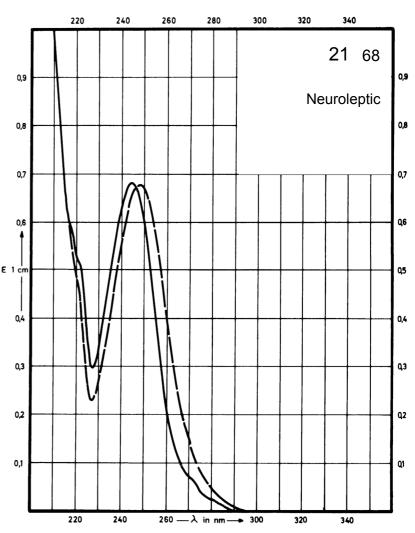

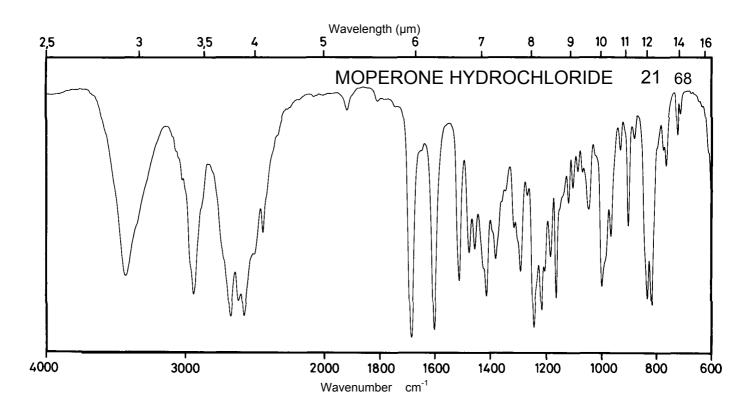

Name MELPERONE HYDROCHLORIDE


 M_r


299.8

Concentration


1.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm	247 nm	247 nm	247 nm
E ^{1%} 1cm	413	416	416	408
ε	12380	12470	12470	12230

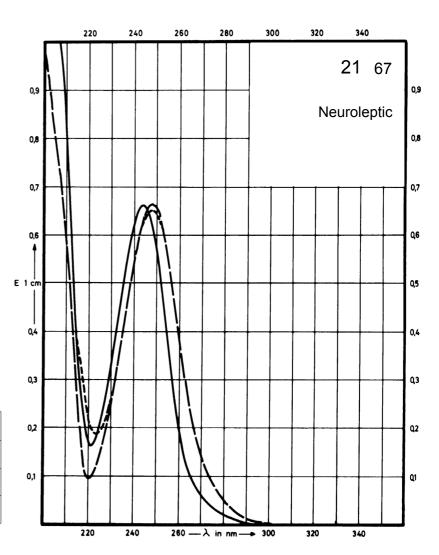

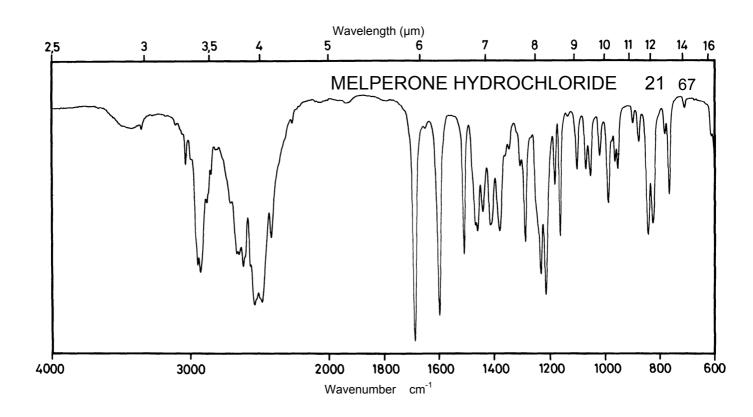




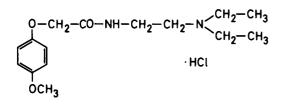


 M_r 316.8


Concentration 2 mg / 100 ml 8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm 224 nm		285 nm 222 nm	285 nm
E ^{1%} 1cm	82 317		73 290	73
3	2600 10030		2310 9190	2310

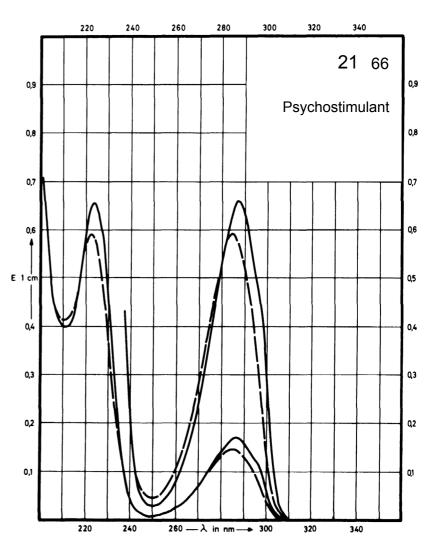

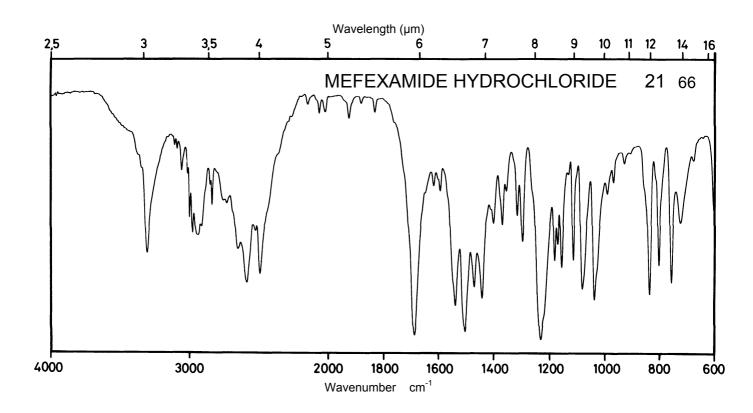





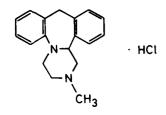


 M_r 300.8

Concentration 1 mg / 100 ml 5 mg / 100 ml

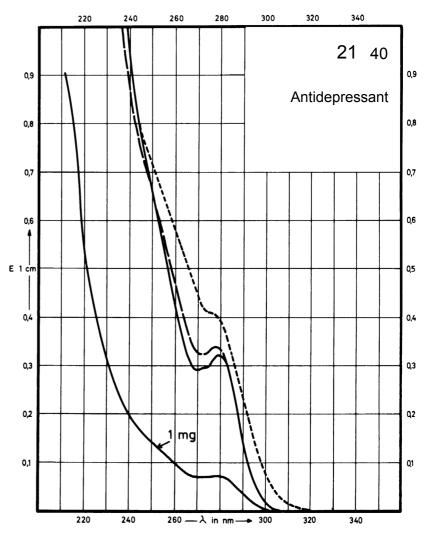

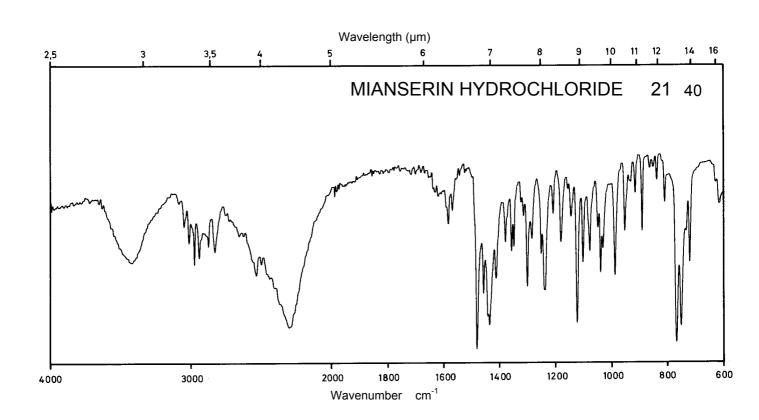




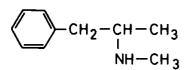


 M_r

Concentration 25 mg / 100 ml

149.2

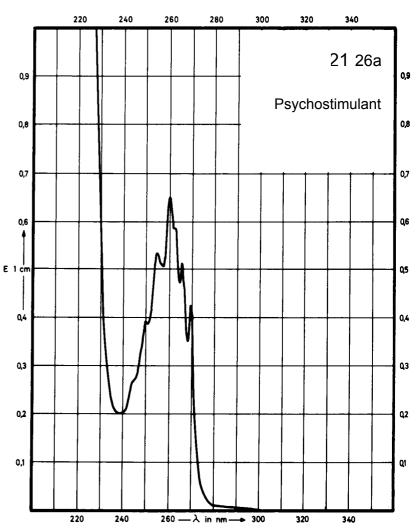

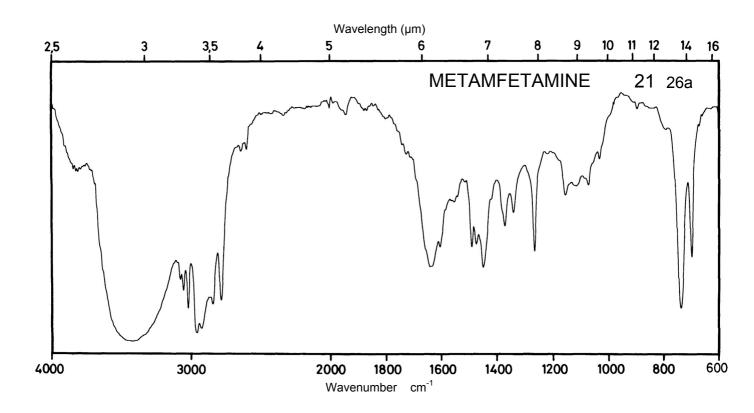


Name METAMFETAMINE HYDROCHLORIDE


 M_{r}

185.7


Concentration 60 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm	263 nm 257 nm 252 nm	263 nm 257 nm 252 nm	
E ^{1%} 1cm	7.5 9.6 7.4	7.6 10.0 8.1	7.6 10.0 8.1	
ε	140 180 138	142 186 150	142 186 150	

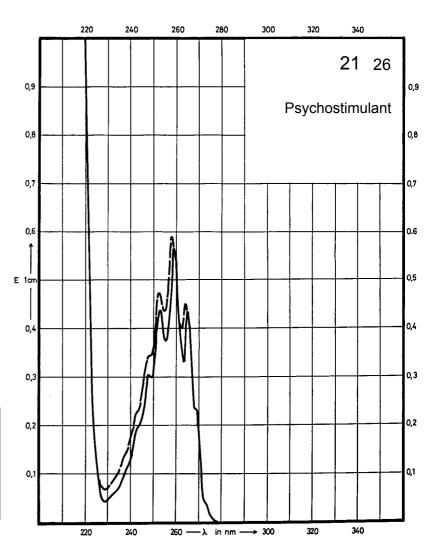

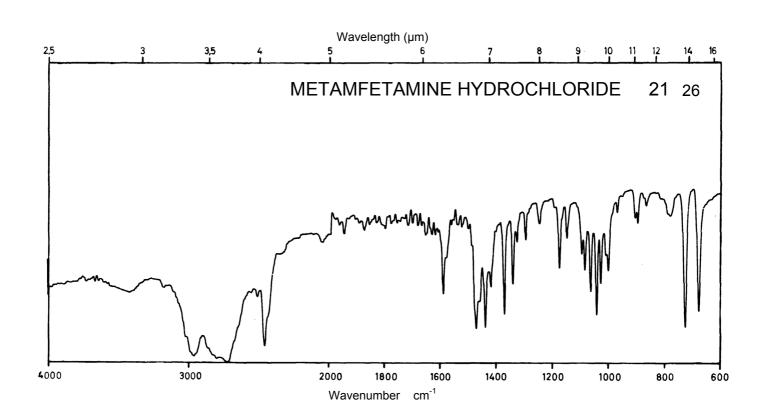

Name METHYLPHENIDATE **HYDROCHLORIDE**

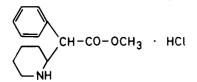


 M_r


233.3

Concentration 50 mg / 100 ml


136

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 264 nm 258 nm 252 nm 5.9 7.4 264 nm 258 nm 264 nm 258 nm 252 nm Maximum of absorption 252 nm 5.9 7.8 5.9 8.0 E^{1%} 1cm 5.8 138 173 6.6 138 187 182 147 ε

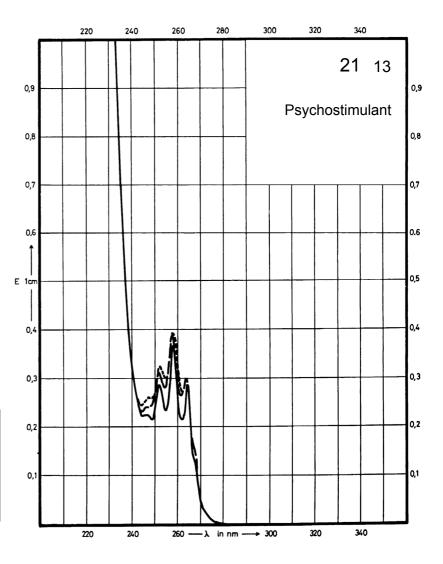

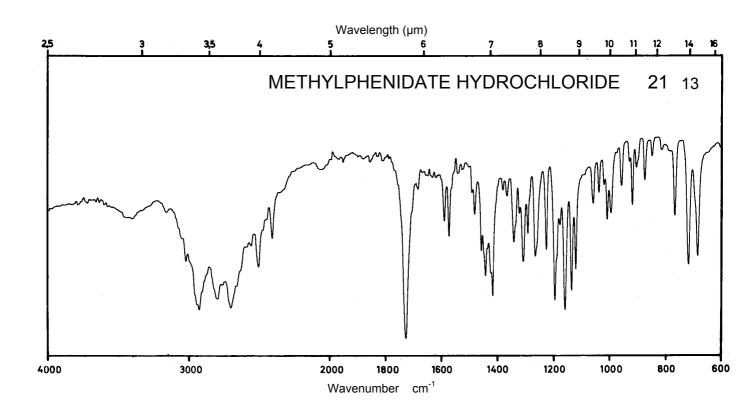


Name MAPROTILINE HYDROCHLORIDE

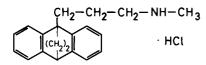


 M_r 313.9


Concentration 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 265 nm	272 nm 265 nm	272 nm 265 nm	
E ^{1%} 1cm	49 40	47 40	47 40	
ε	1540 1250	1480 1250	1480 1250	

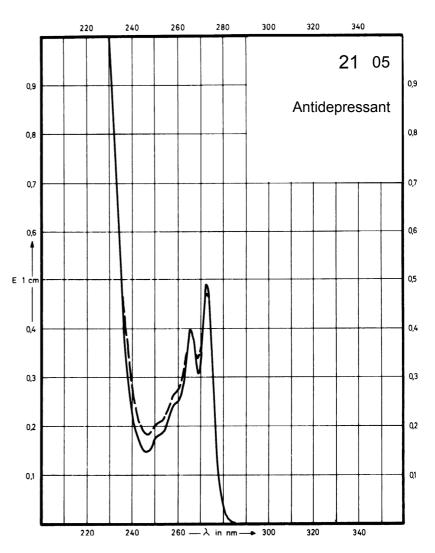

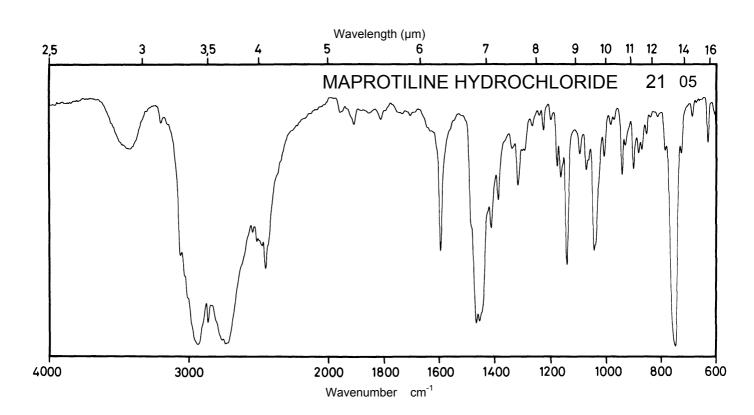




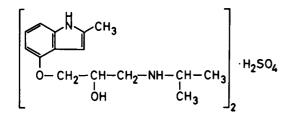


 M_r 622.8


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm 267 nm	267 nm	267 nm	267 nm
E ^{1%} 1cm	149 328	303	303	307
ε	9300 20400	18900	18900	19100

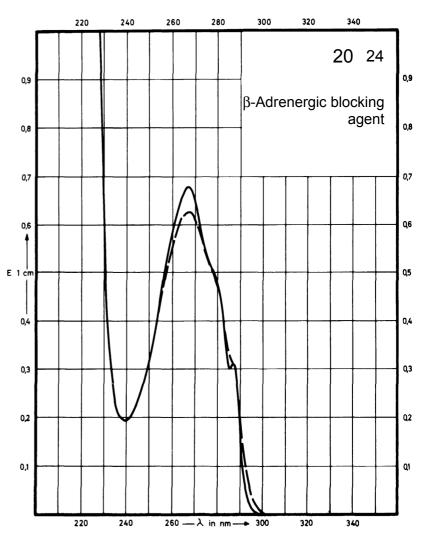

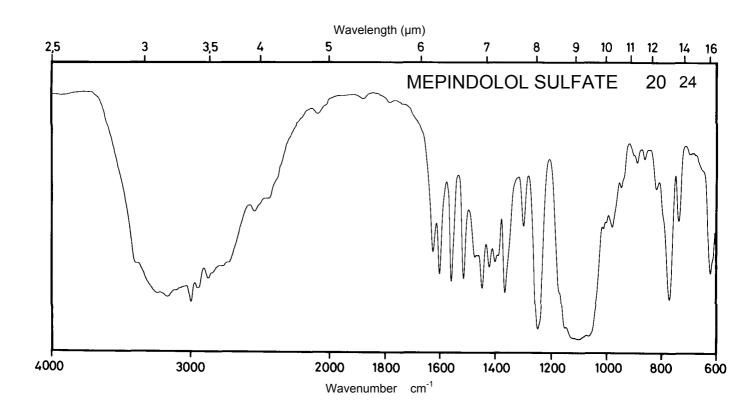




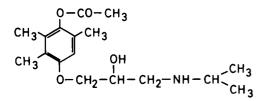


 M_{r} 309.4

Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	279 nm		278 nm	Decom- position observed
E ^{1%} 1cm	53		44	
ε	1640		1360	

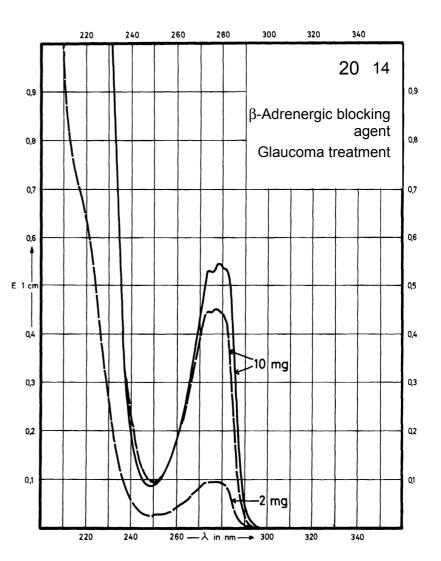

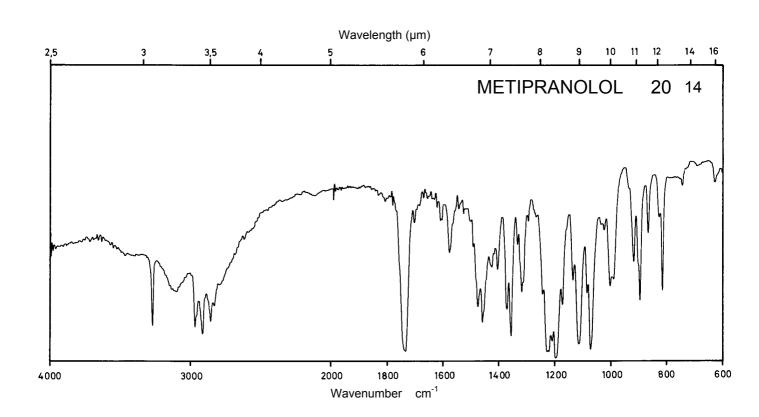

Name METOPROLOL TARTRATE


$$\begin{array}{c} \text{CH}_3 \\ \text{O-CH}_2\text{-CH-CH}_2\text{-NH-CH} \\ \text{OH} \\ \text{CH}_3 \\ \text{CH}_2\text{-CH}_2\text{-O-CH}_3 \end{array}$$


 M_r

417.5


Concentration 2.5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 276 nm 224 nm		274 nm 222 nm	275 nm
E ^{1%} 1cm	38.0 45.2 304		40.9 278	41.3
3	1590 1890 12690		1710 11610	1730

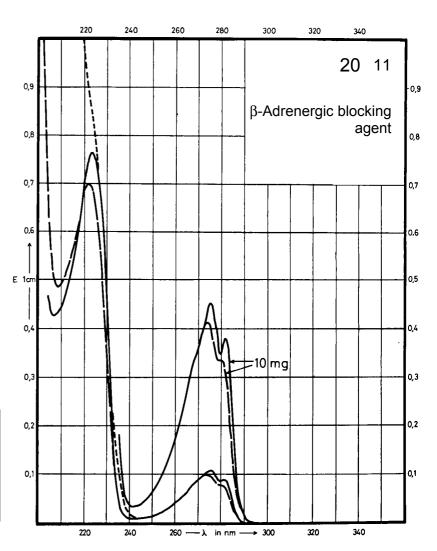

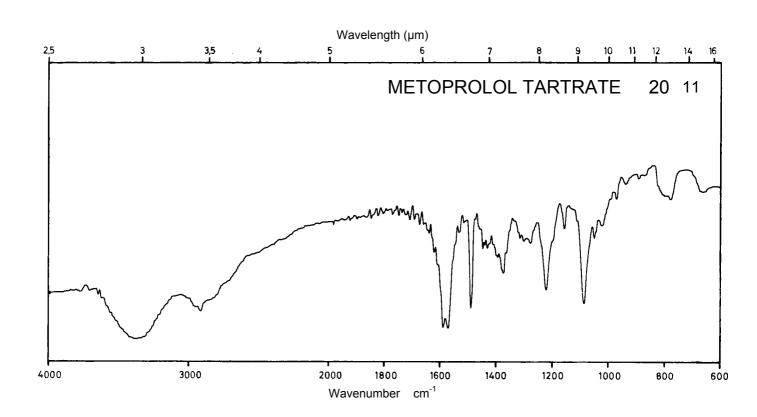




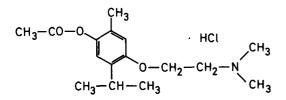


 M_r 315.8


Concentration 1 mg / 100 ml 5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm		275 nm	301 nm
E ^{1%} 1cm	70		70	129
ε	2210		2210	4070

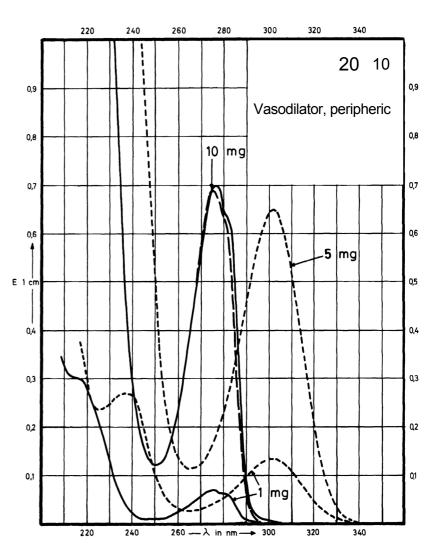

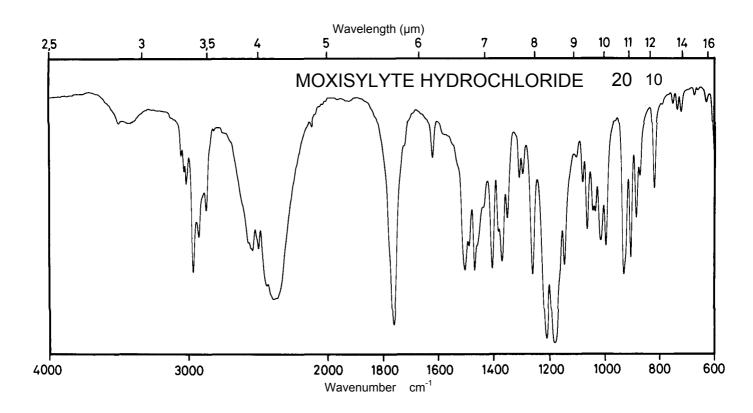



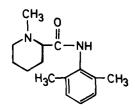



246.3

Concentration


52 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 263 nm		271 nm 262 nm	271 nm 263 nm
E ^{1%} 1cm	10.3 14.2		15.7 19.0	13.4 16.8
ε	250 350		390 470	330 410

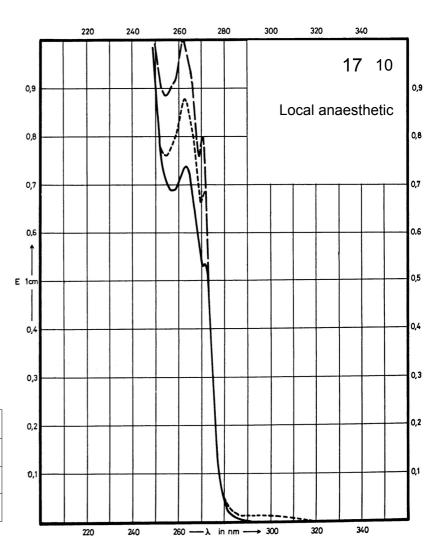

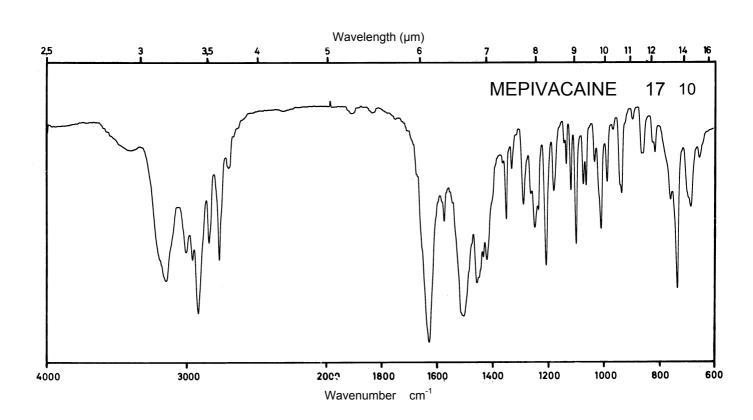



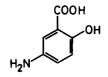



 M_{r}

153.1

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			302 nm 231 nm	331 nm
E ^{1%} 1cm			229 445	222
ε			3500 6800	3400

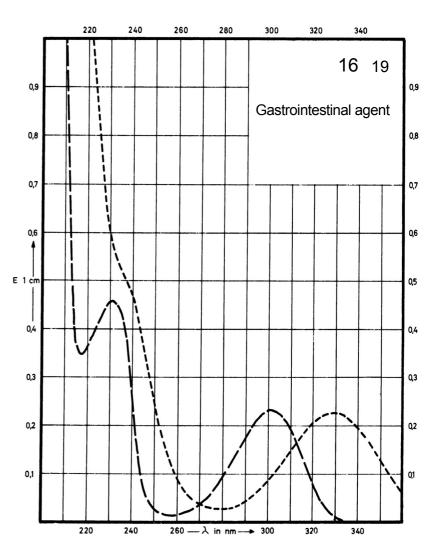

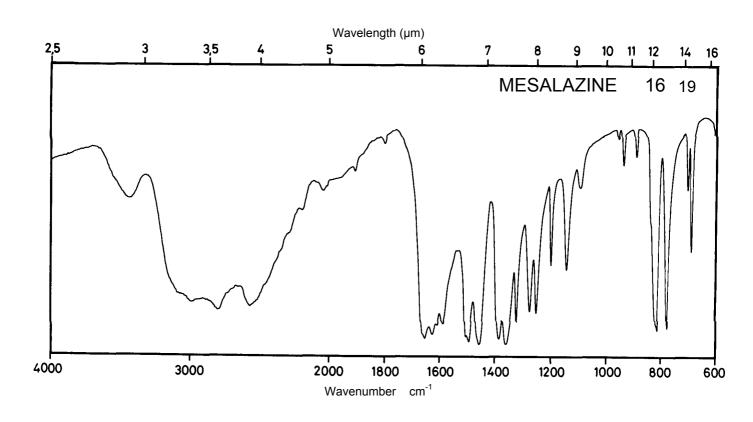


Name METOCLOPRAMIDE HYDROCHLORIDE


$$\begin{array}{c} \text{OC-NH-CH}_2\text{-CH}_2\text{-N} \\ \text{OCH}_3 \\ \text{CI} \\ \text{NH}_2 \end{array} \\ \begin{array}{c} \text{CH}_2\text{-CH}_3 \\ \text{CH}_2\text{-CH}_3 \\ \end{array}$$

 M_r 336.3


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm 277 nm 213 nm		309 nm 273 nm 212 nm	308 nm 271 nm
E ^{1%} 1cm	441 487 810		315 345 850	333 420
ε	14830 16380 27240		10590 11600 28600	11200 14120

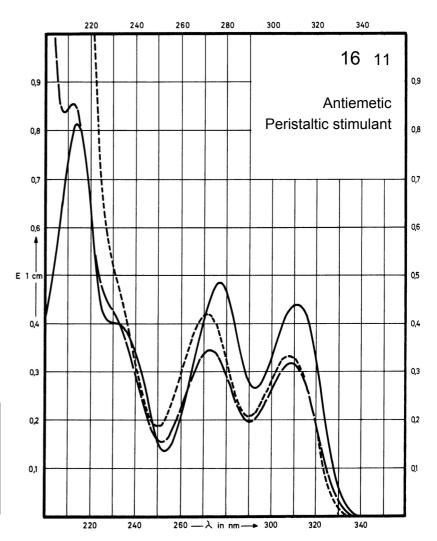

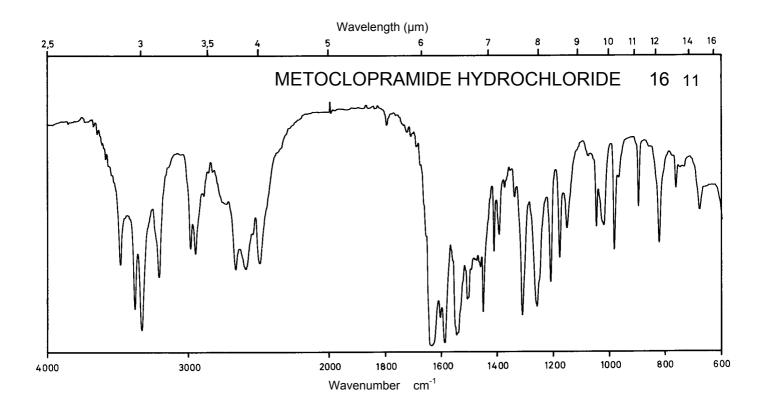




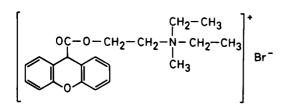


420.4


Concentration 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm 247 nm	282 nm	282 nm	283 nm 250 nm
E ^{1%} 1cm	67 133	64	64	62 133
3	2830 5600	2700	2700	2620 5600

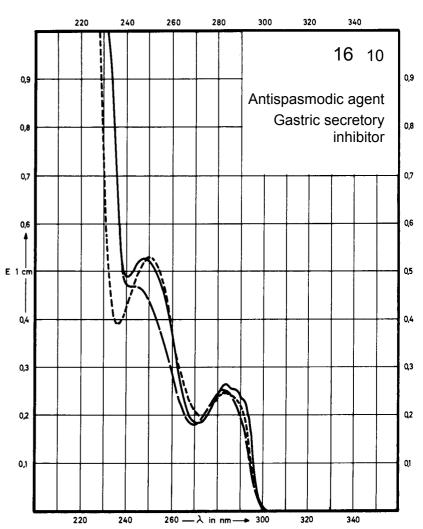

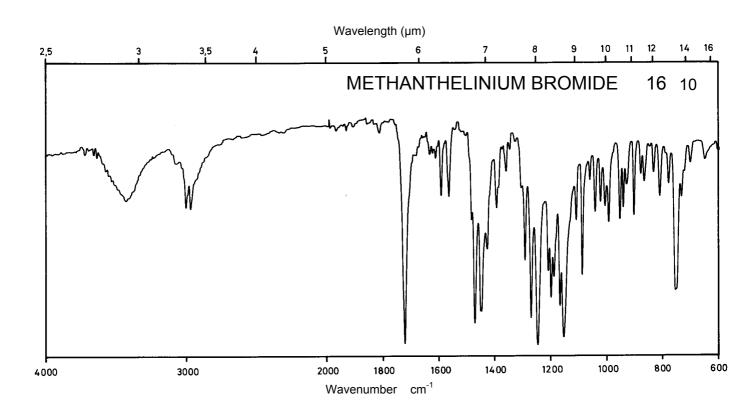



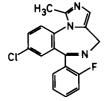



 M_r

Concentration 0.5 mg / 100 ml


325.8

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol	-	:		
Maximum of absorption				
E ^{1%} 1cm				
3				

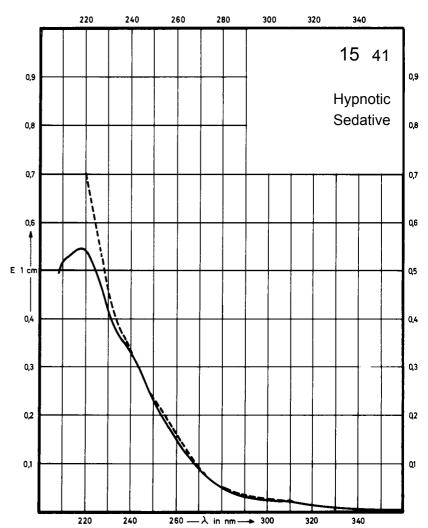

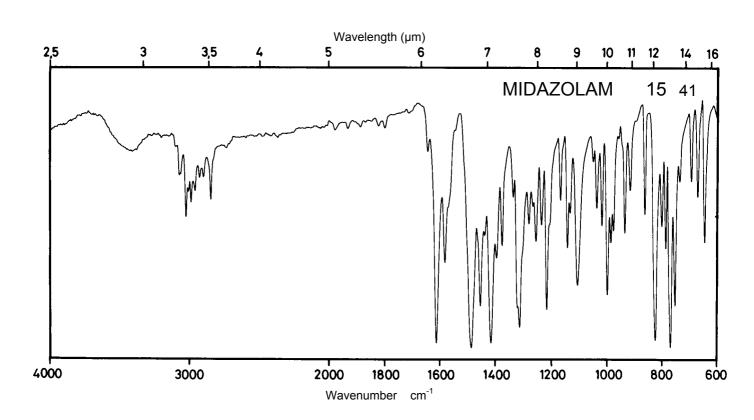

Name METHOHEXITAL

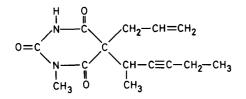



 M_{r}

262.3


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	218 nm		222 nm	247 nm
E ^{1%} 1cm	238		253	290
ε	6200		6600	7600

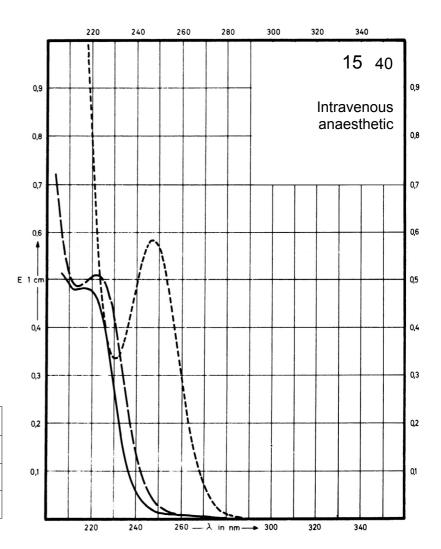

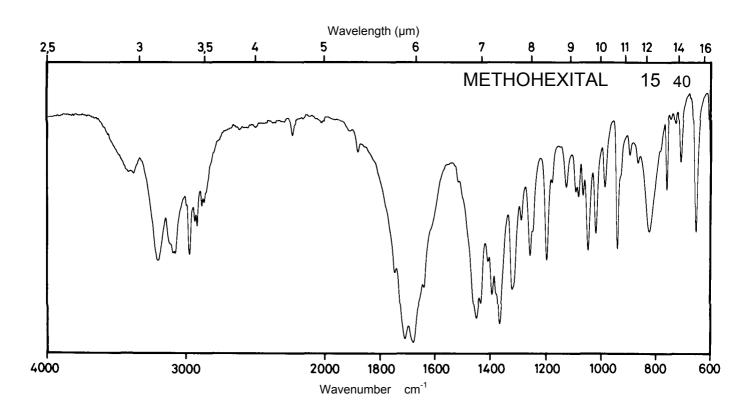


Name MESUXIMIDE

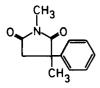


 M_r 203.2


Concentration 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm 248 nm			
E ^{1%} 1cm	18.4 22.5			
ε	373 457			

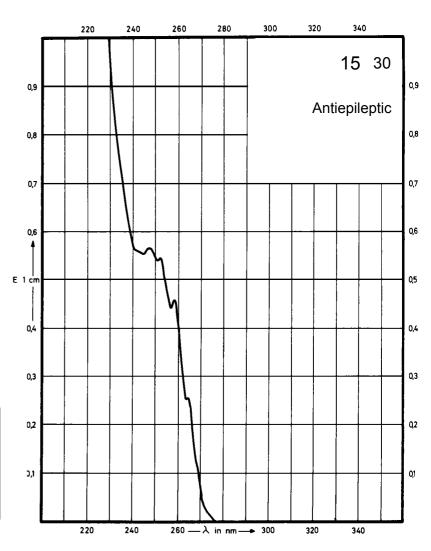

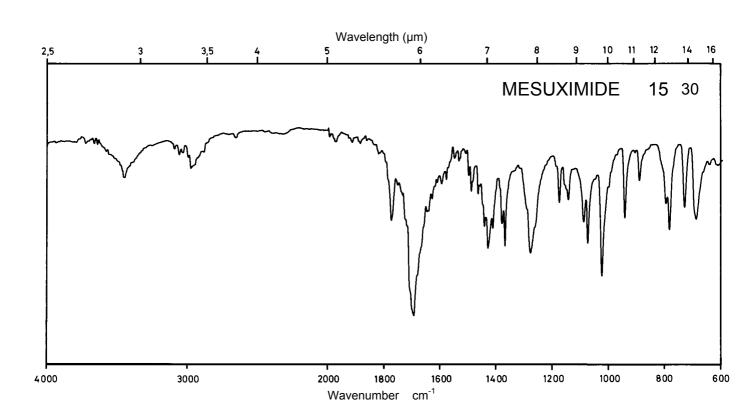




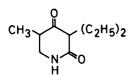


183.3


Concentration 50 mg / 100 ml 200 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	304 nm		300 nm	375 nm 286 nm
E 1% 1cm	4.30		3.84	3.68 2.40
ε	79		70	67 44

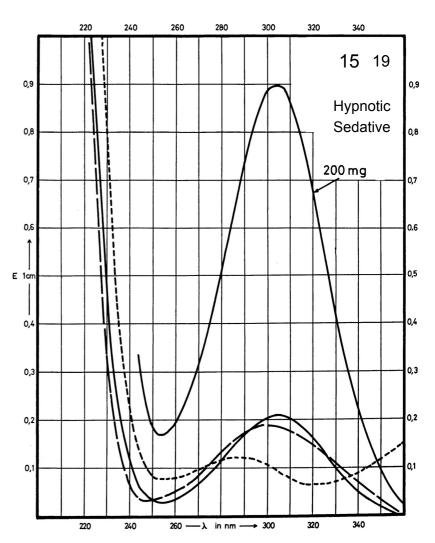

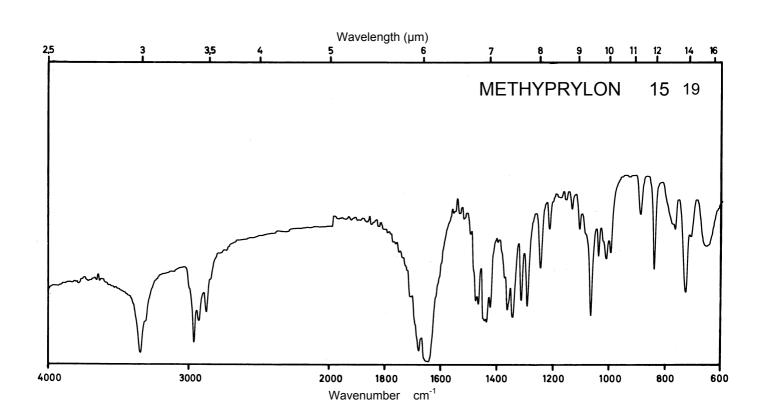

Name METHYLPHENO-BARBITAL

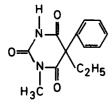



 M_{r}

246.3


Concentration 2 mg / 100 ml 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				245 nm
E ^{1%} 1cm				343
3				8450

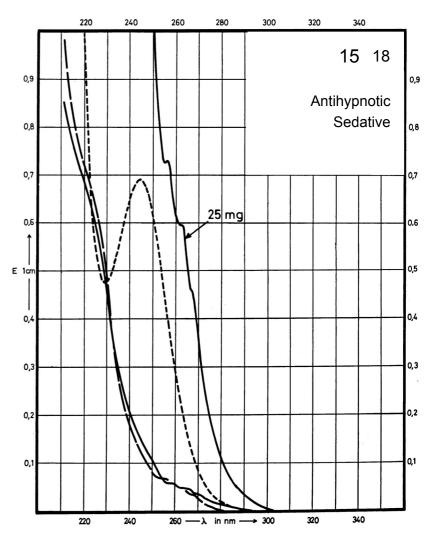

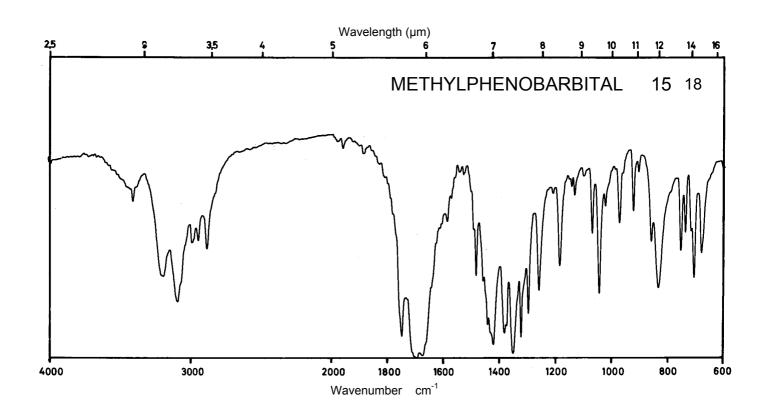

Name METHAQUALONE

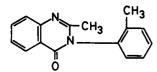


 M_{r}


250.3

Concentration


2.3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	304 nm 264 nm		269 nm	305 nm 264 nm
E ^{1%} 1cm	149 392		319	154 347
ε	3730 9810		7980	3850 8690

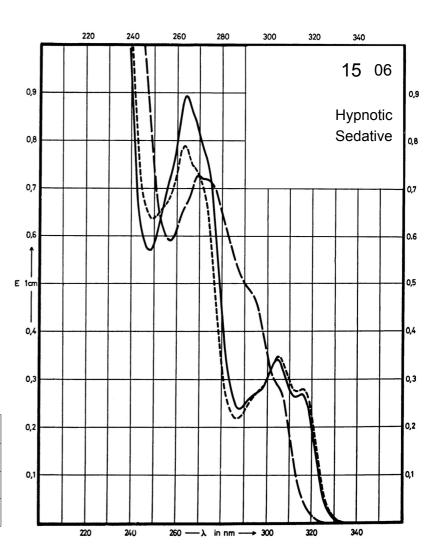

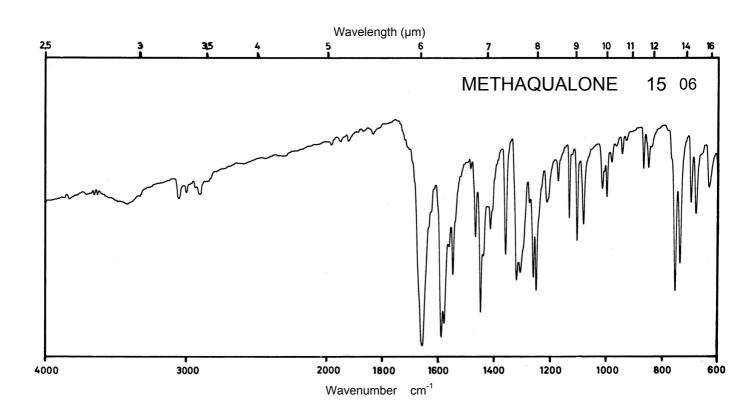




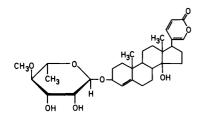


544.7


Concentration 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	299 nm		300 nm	
E ^{1%} 1cm	107		119	
ε	5800		6500	

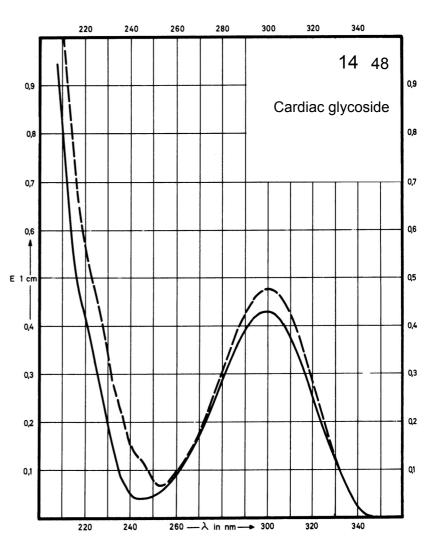

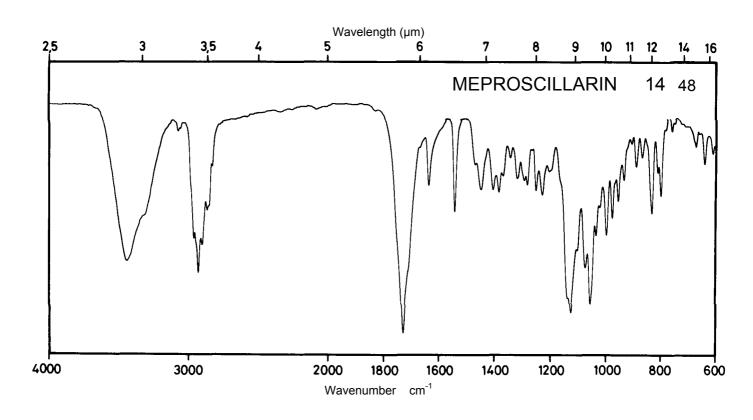




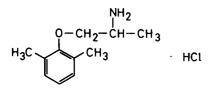


 $\,M_r\,$

215.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm	262 nm	262 nm	
E ^{1%} 1cm	11.1	11.8	11.8	
ε	240	254	254	

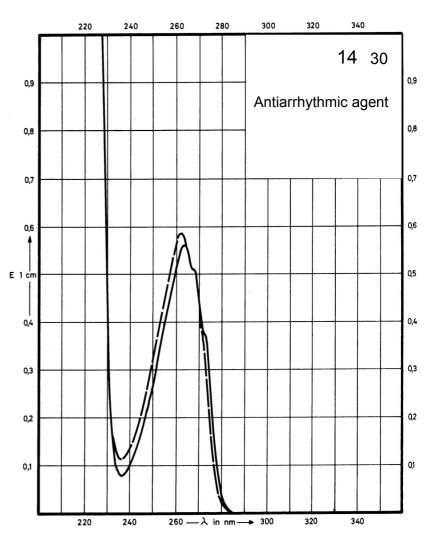

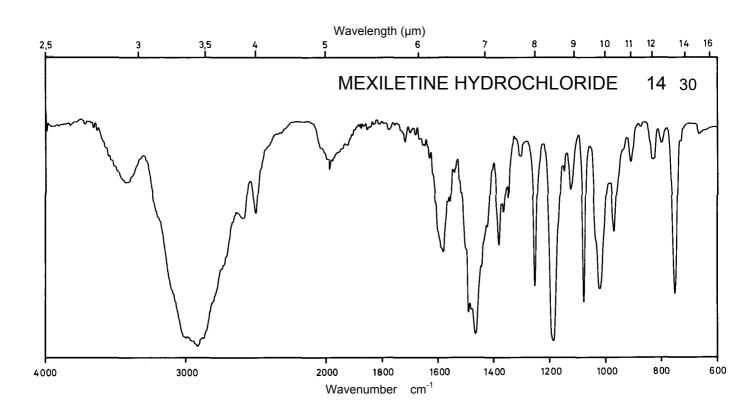



Name MOLSIDOMINE

 M_{r}


242.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 228 nm		286 nm	310 nm
E ^{1%} 1cm	580 580		684	625
ε	14050 14050		16560	15140

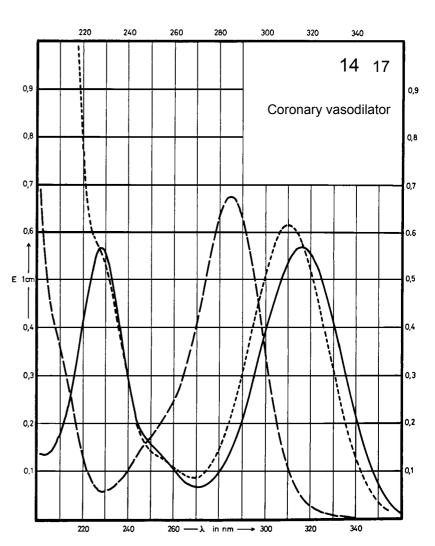

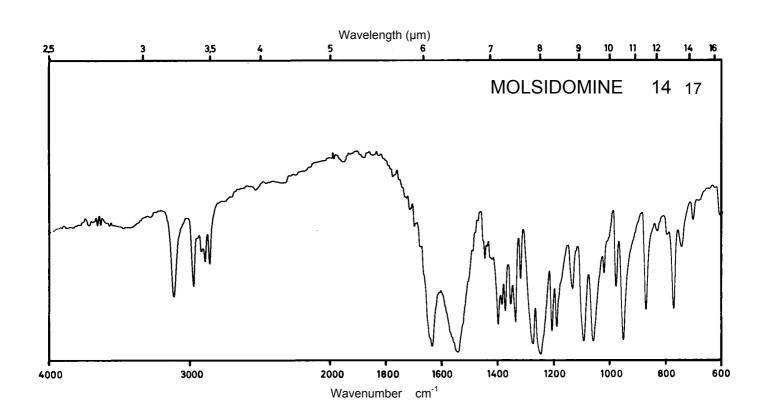



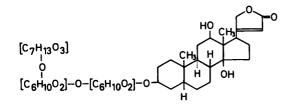



795.0

Concentration


2.7 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	221 nm		218 nm	
E ^{1%} 1cm	179		222	
ε	14230		17650	

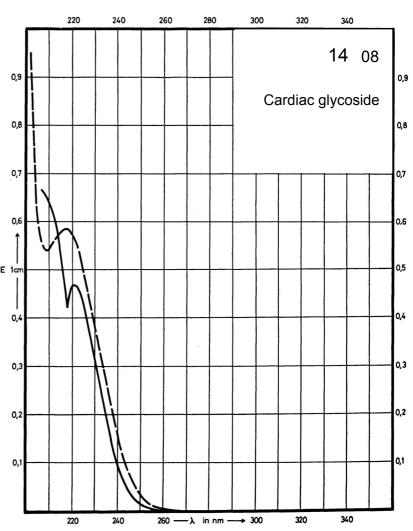

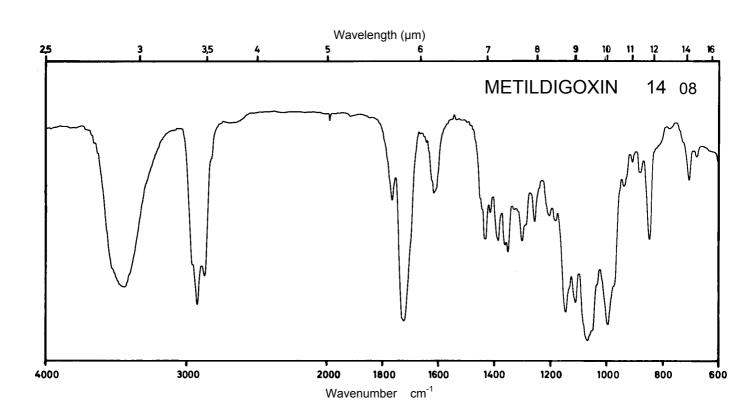



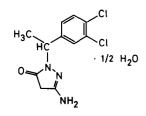


 $\,M_r\,$


272.1 (anhydrous)

Concentration


2.4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm		271 nm	270 nm 280 nm
E ^{1%} 1cm	154		139	37 28
3	4175		3770	996 760

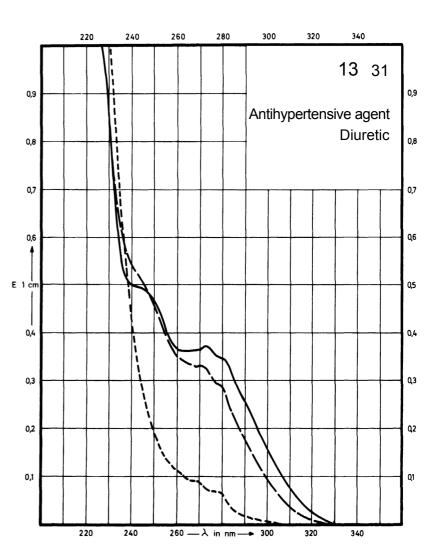

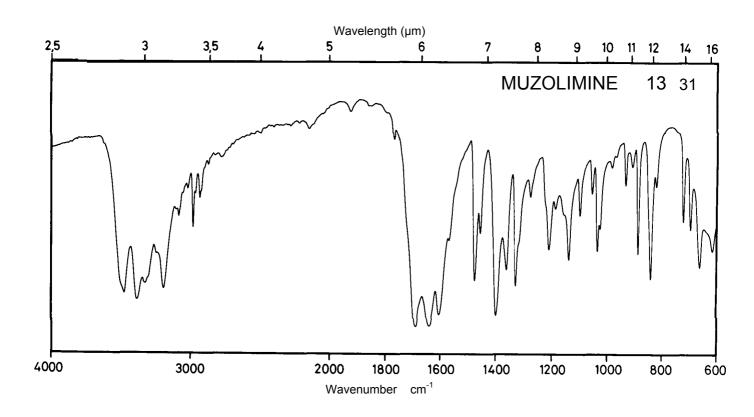





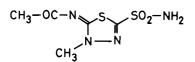


236.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 254 nm	290 nm 252 nm	290 nm 255 nm	288 nm
E ^{1%} 1cm	523 197	521 202	416 242	579
ε	12360 4655	12310 4770	9830 5720	13680

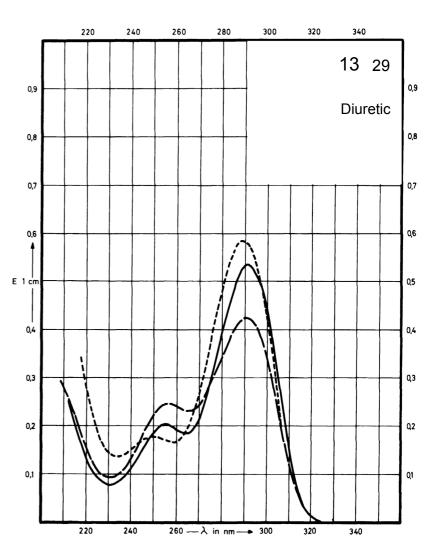

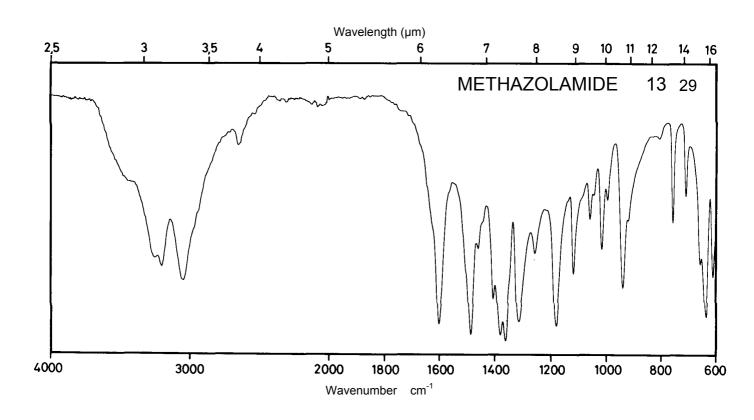



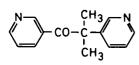



226.3

Concentration

2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 230 nm		260 nm	264 nm
E ^{1%} 1cm	265 350		240	143
ε	5990 7930		5420	3230

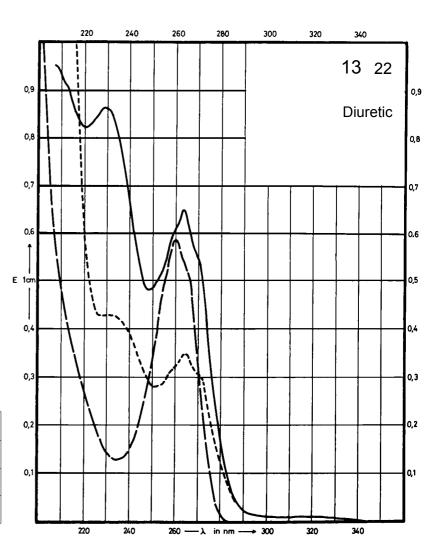

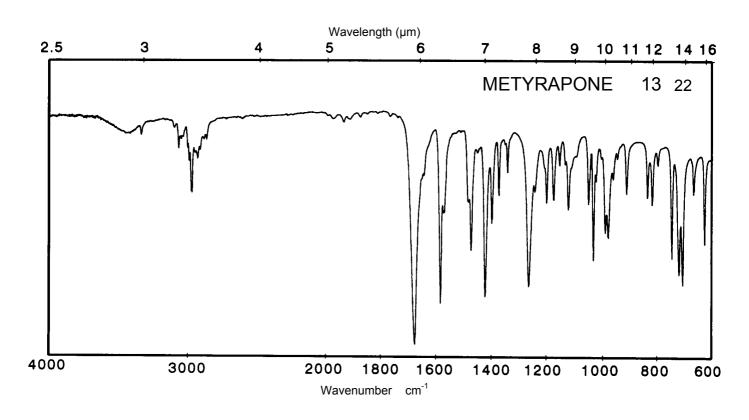

Name METOLAZONE


H₂NO₂S O CH₃


 M_{r}

365.8

Concentration 0.5 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	343 nm 271 nm 236 nm		344 nm 271 nm 237 nm	347 nm 236 nm
E 1% 1cm	89 317 1530		81 298 1470	79 1500
ε	3250 11600 55980		2970 10890 53720	2880 54890

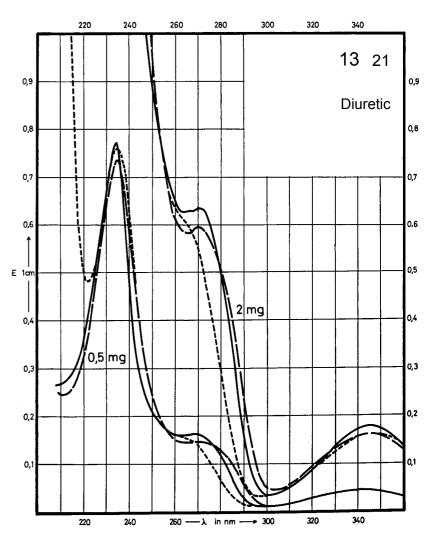

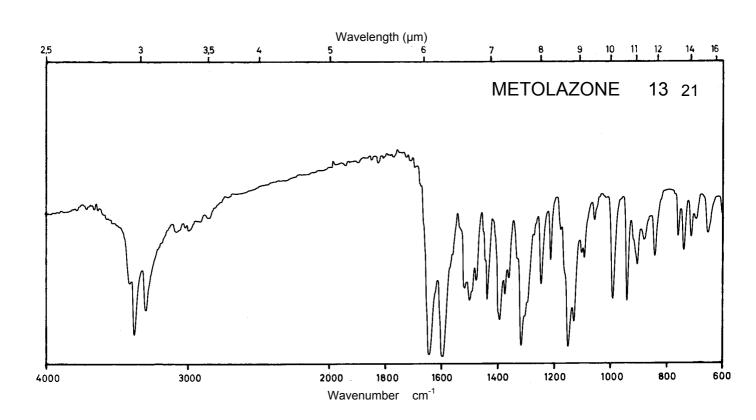





382.9


Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm 276 nm		284 nm 276 nm	
E ^{1%} 1cm	31 44		38 50	
ε	1190 1680		1460 1910	

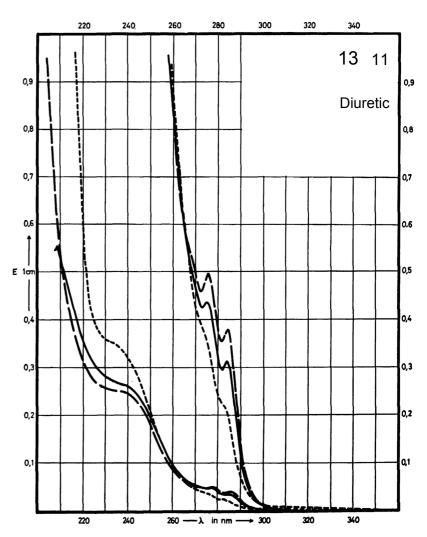

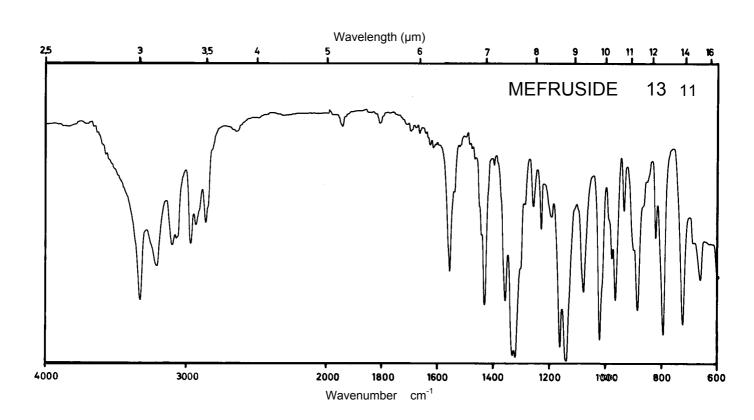




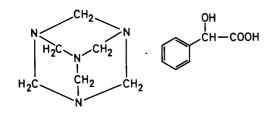


 M_r


292.3

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	264 nm		262 nm	
	258 nm		256 nm	
absorption	252 nm		251 nm	
1%	5.09		5.90	
E ^{1%} 1cm	6.72		7.30	
1cm	5.50		5.80	
	150		170	
3	195		215	
	160	l	170	

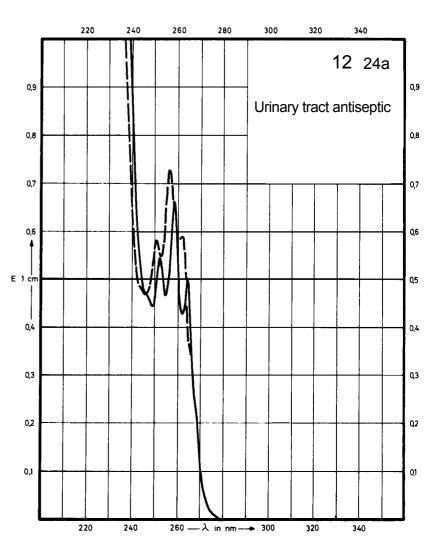

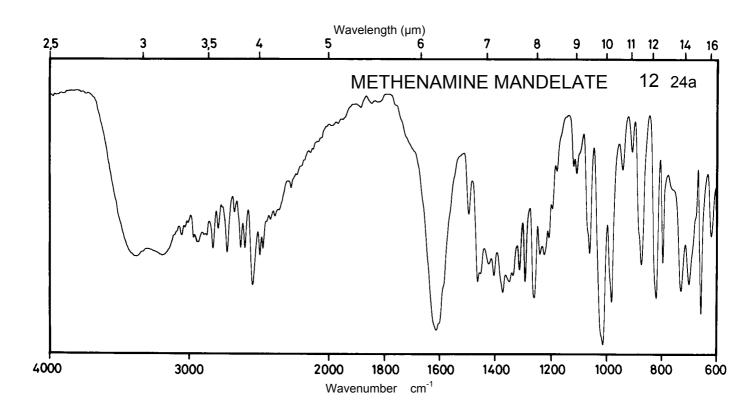





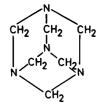


 M_r 140.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

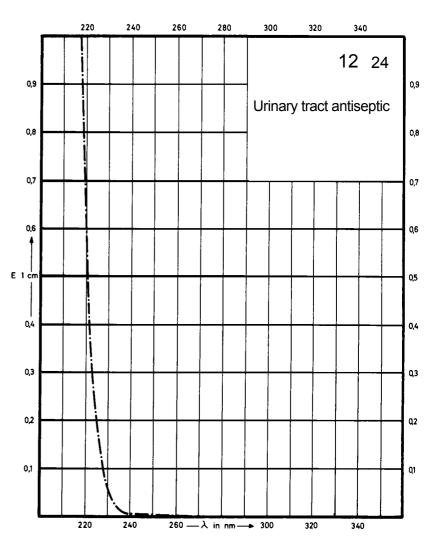

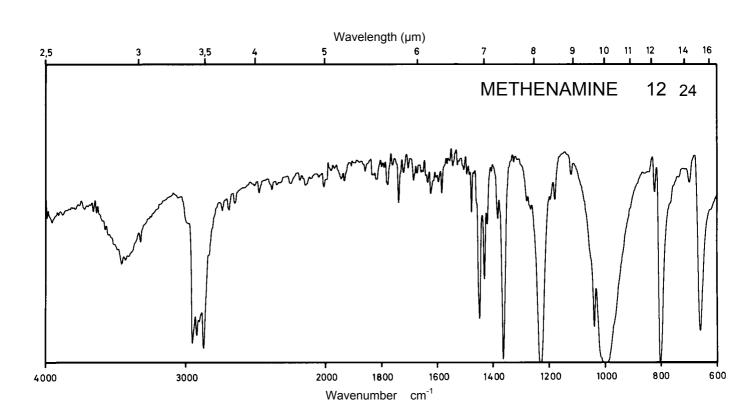



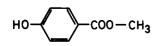



152.1

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	256 nm		255 nm	295 nm
E ^{1%} _{1cm}	1063		1006	1447
ε	16170		15300	22010

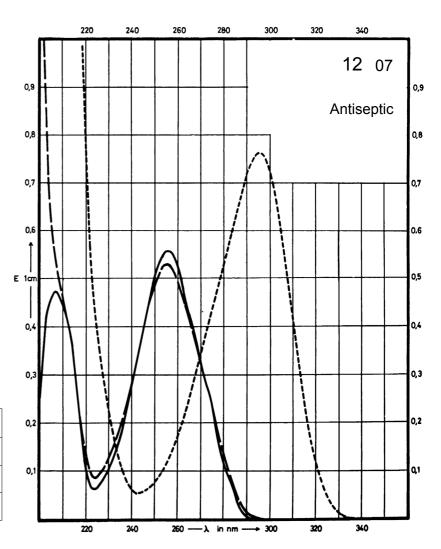

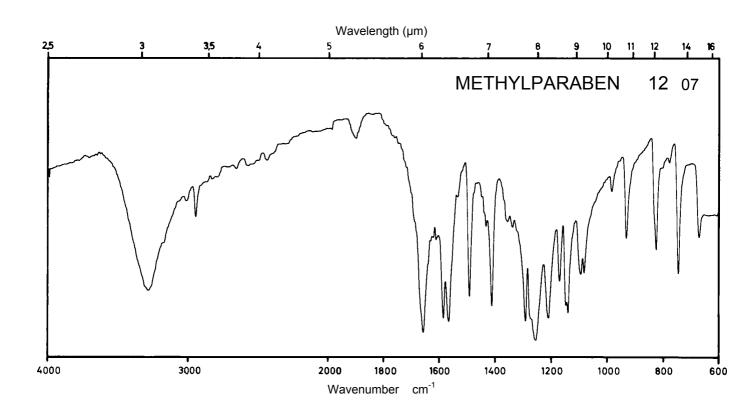



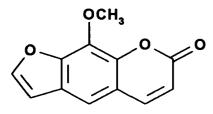



216.2


Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	299 nm 248 nm 218 nm	302 nm 247 nm 217 nm	303 nm 247 nm 217 nm	345 nm 250 nm
E ^{1%} 1cm	567 1121 1158	558 1003 1105	559 1000 1105	245 984
ε	12200 24200 25000	12100 21700 23900	12100 21600 23900	5300 21300

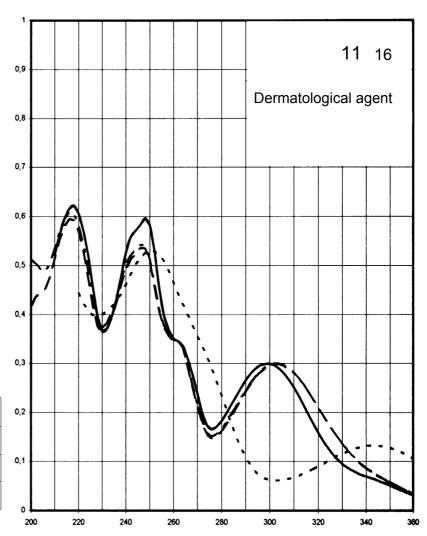

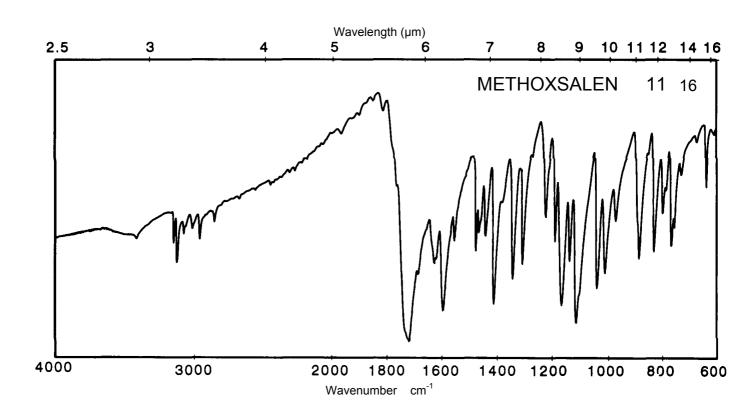

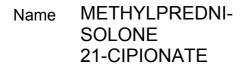


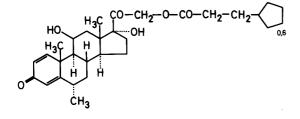


 M_r


498.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm			
E ^{1%} 1cm	302			
ε	15000			

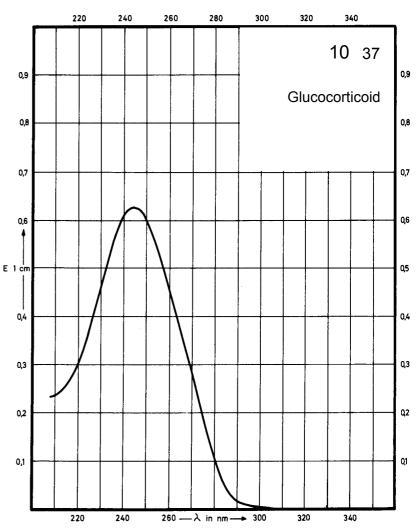

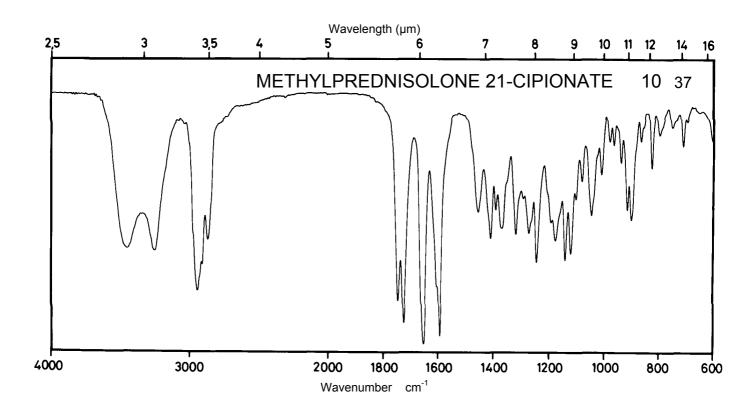



Name METHYLPREDNI-SOLONE 21-(SODIUM SUCCINATE)

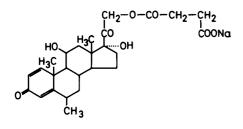


 M_r 496.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	244 nm	248 nm		
E ^{1%} 1cm	296	291		
3	14700	14440		

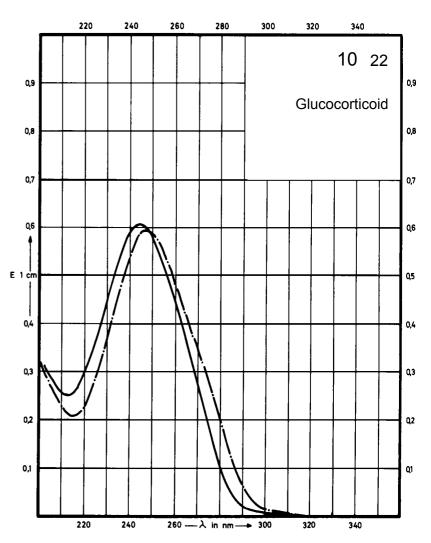

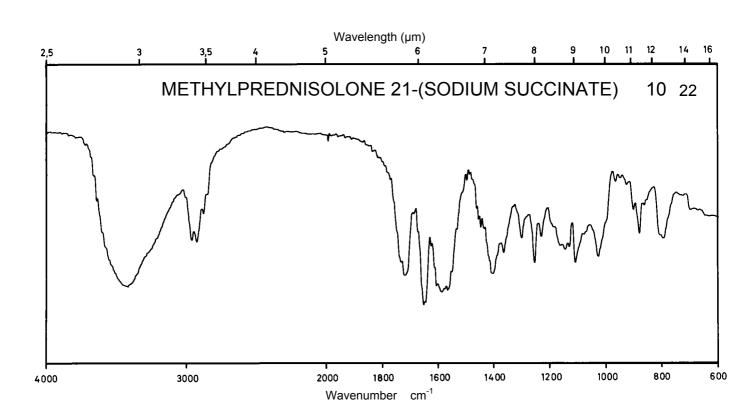



Name METHYLPREDNI-SOLONE 21-ACETATE

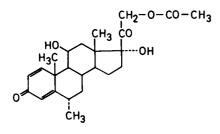


 M_r 416.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm	247 nm		
E ^{1%} 1cm	360	360		
3	14990	14990		

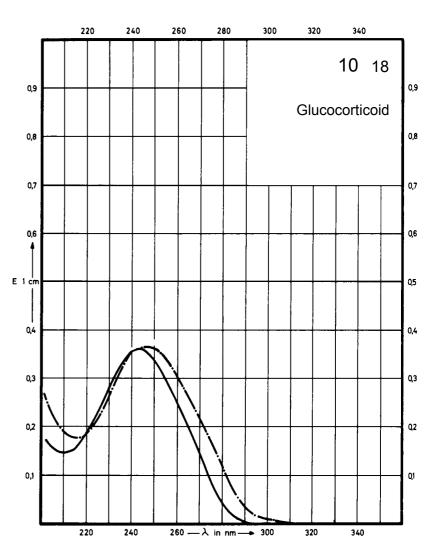

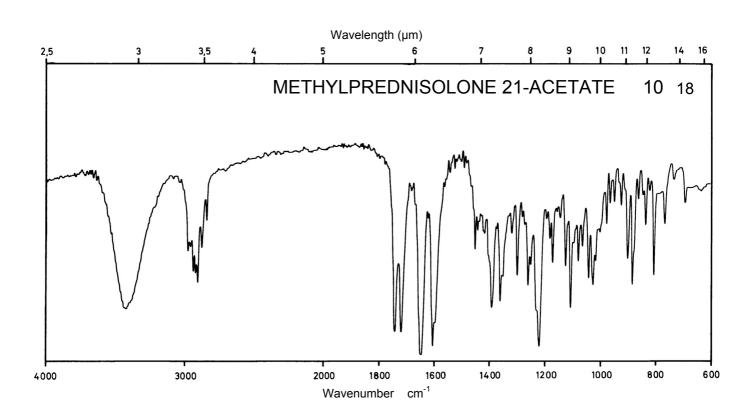


Name METANDIENONE

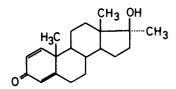


 M_{r}


284.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm		250 nm	250 nm
E ^{1%} 1cm	520		520	520
ε	14800		14800	14800

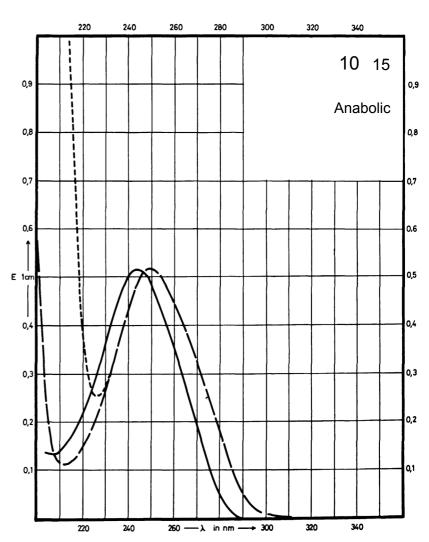

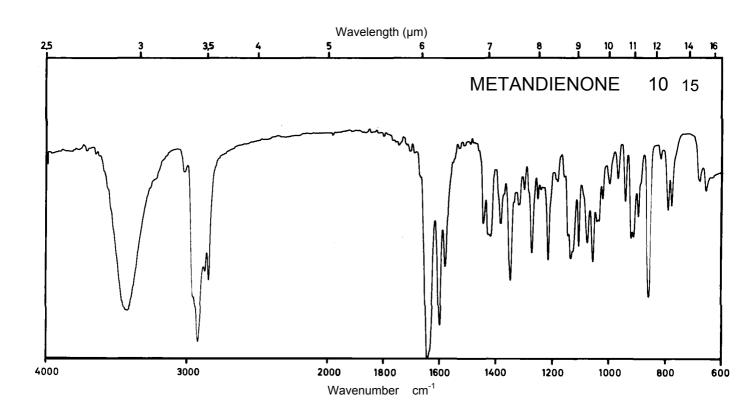


Name METHYLPREDNI-SOLONE

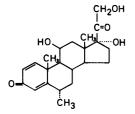


 M_{r}


374.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	242 nm		246 nm	246 nm
E ^{1%} 1cm	404		399	399
ε	15130		14940	14940

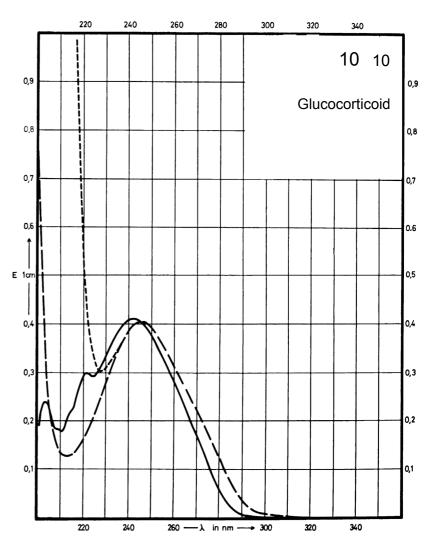

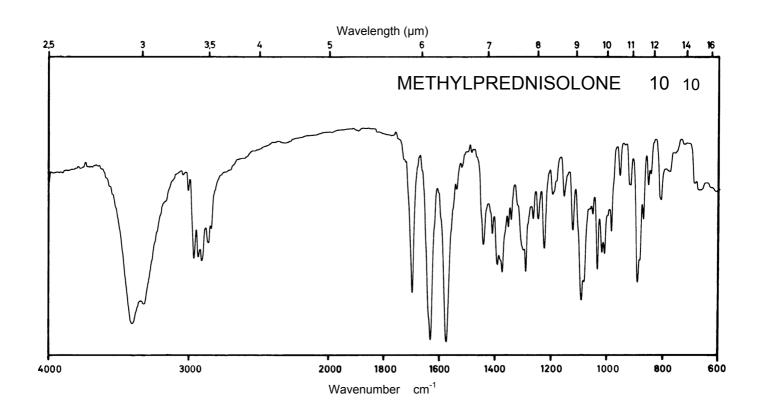





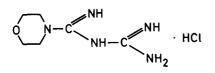


 M_r 207.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm			233 nm
E ^{1%} 1cm	865			675
ε	18000			14000

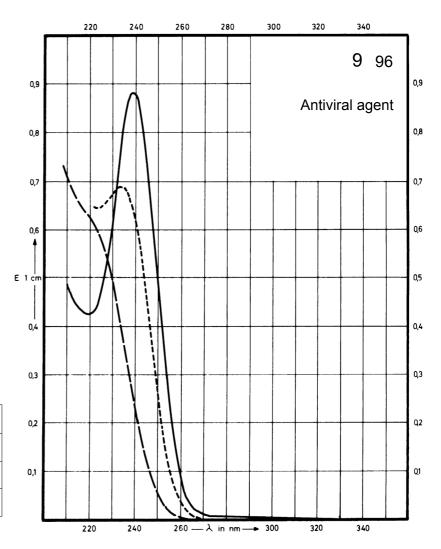

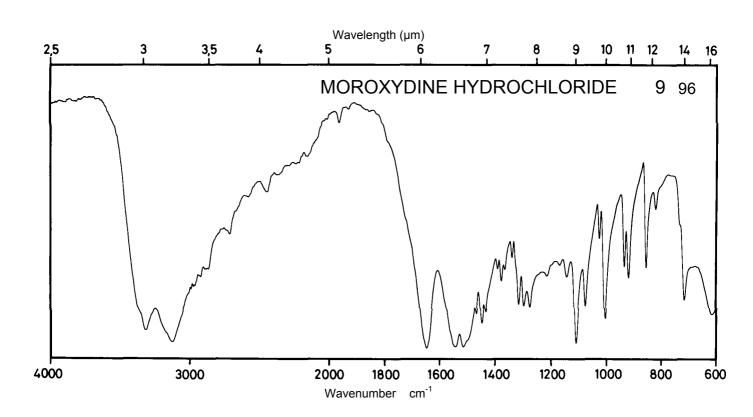





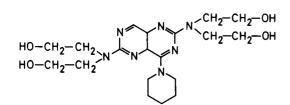


421.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm	285 nm	278 nm 235 nm	287 nm
E ^{1%} 1cm	910	787	804 608	860
3	38400	33200	33900 25600	36200

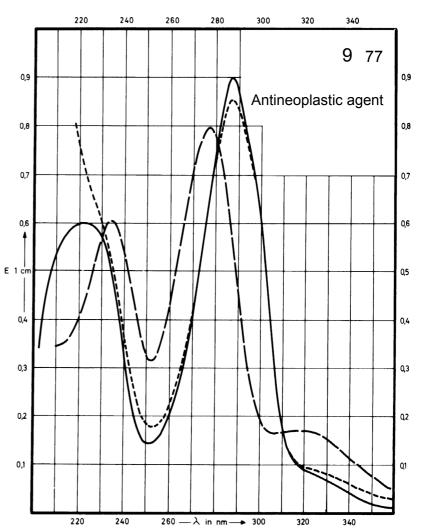

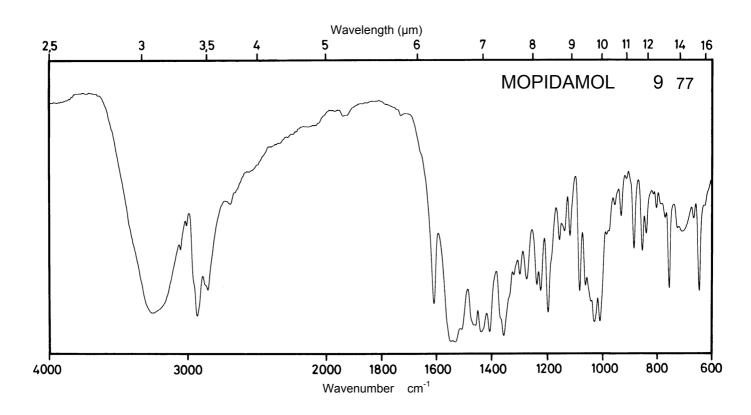





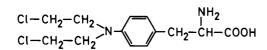

305.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	302 nm 261 nm	302 nm 261 nm	261 nm	303 nm 259 nm
E ^{1%} 1cm	64 707	59 575	146	56 527
ε	1950 21600	1800 17550	4450	1700 16100

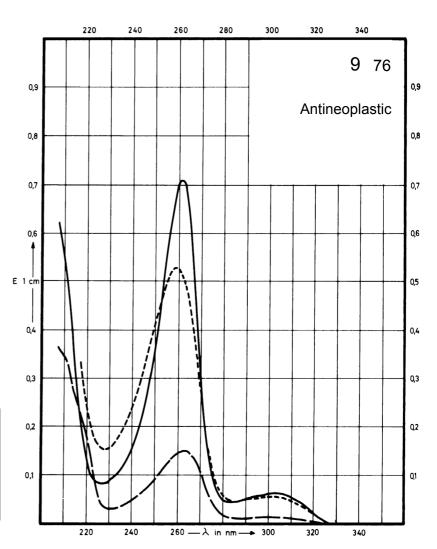

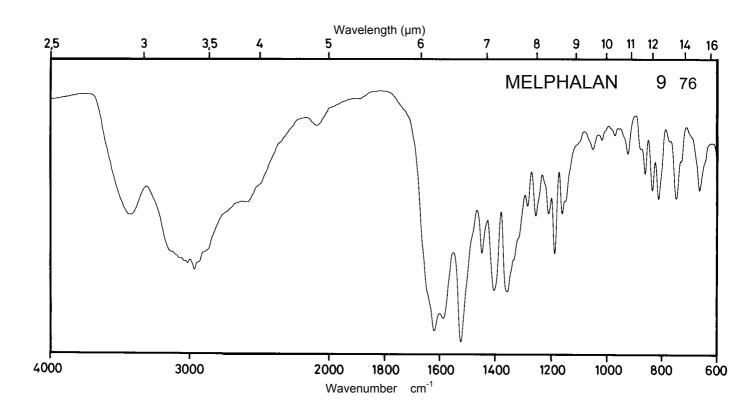

Name MERCAPTOPURINE

SH H


 M_{r}


152.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	328 nm		326 nm	310 nm
E ^{1%} 1cm	1110		1140	1170
3	16890		17350	17810

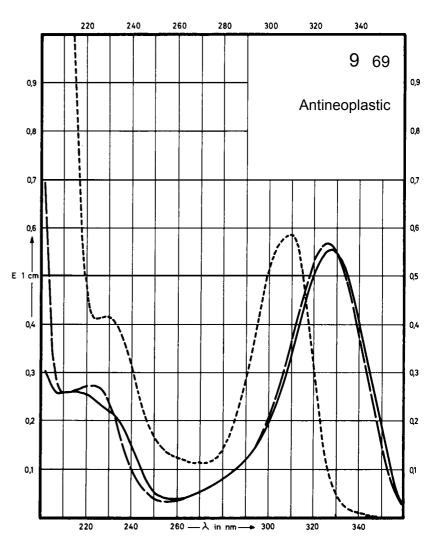

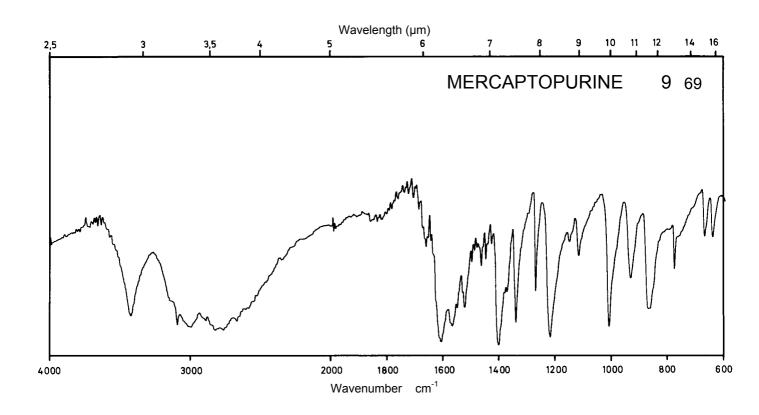





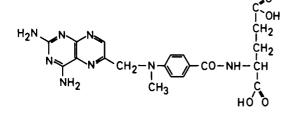

 M_r 454.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	299 nm		306 nm 243 nm	372 nm 303 nm 258 nm
E ^{1%} 1cm	536		439 374	163 498 503
3	24370		19940 17000	7420 22630 22850

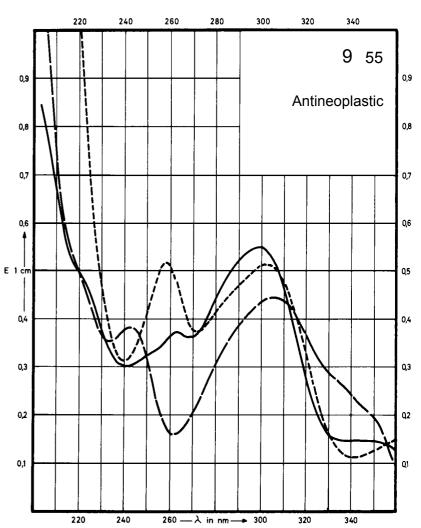

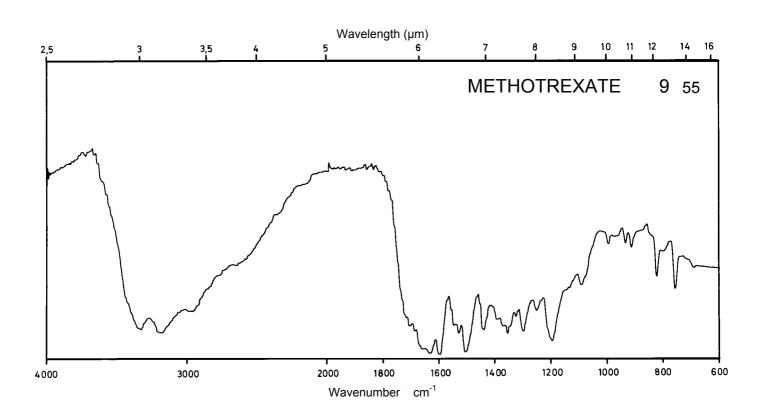





 M_r 479.1


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm 272 nm			
E ^{1%} 1cm	10.9 13.3			
ε	520 640			

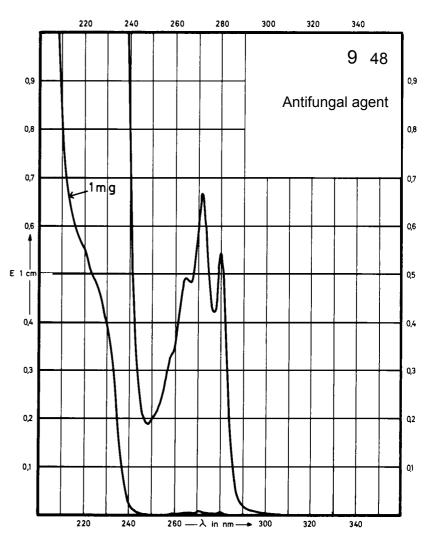

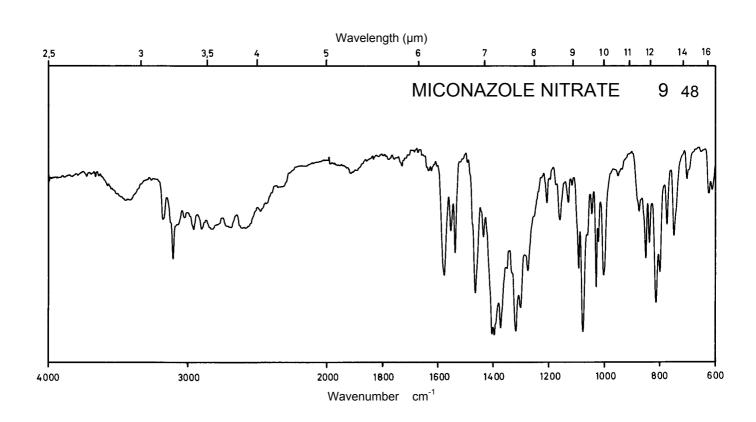


Name METRONIDAZOLE

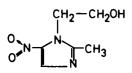


 M_{r}


171.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm		277 nm	319 nm
E ^{1%} 1cm	506		374	520
ε	8660		6400	8900

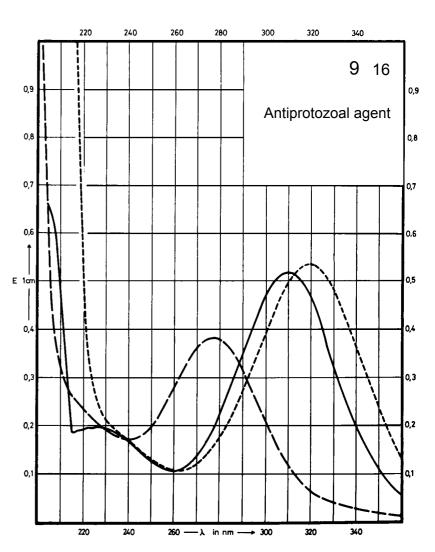

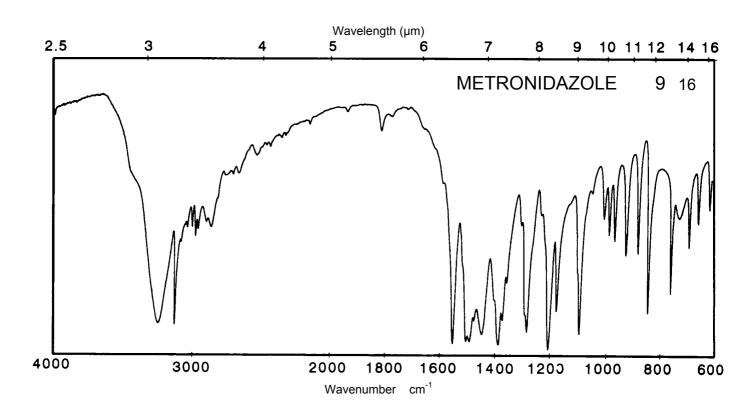


Name MECLOFENAMIC ACID SODIUM SALT

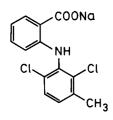


 $\,M_r\,$

318.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	319 nm 289 nm			318 nm 278 nm
E ^{1%} 1cm	200 225			144 190
ε	6400 7200			4900 6000

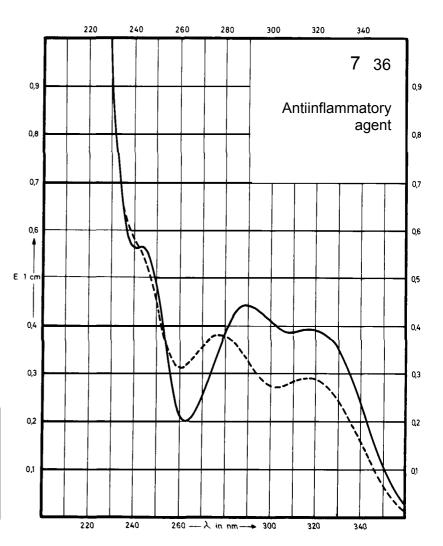

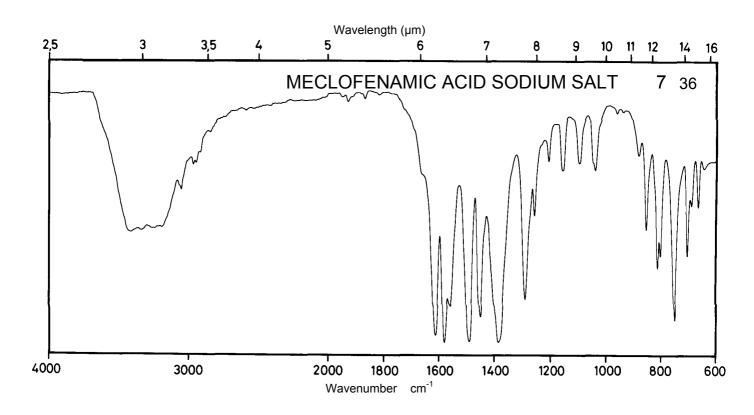



Name MOFEBUTAZONE

 $\,M_r\,$

232.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 240 nm		268 nm 235 nm	279 nm 246 nm
E ^{1%} 1cm	284 495		189 456	326 418
ε	6600 11500		4390 10600	7580 9700

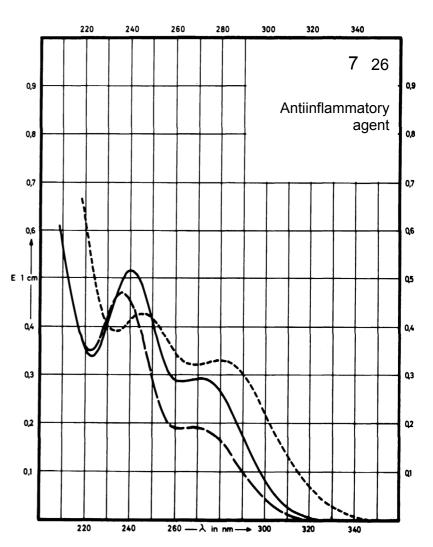

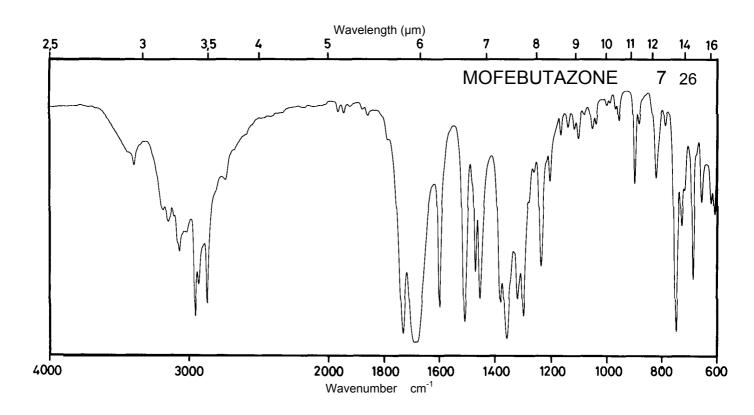

Name MEFENAMIC ACID

COOH

CH₃ CH₃ CH₃

 M_{r}


241.3


Concentration 1

1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI*	0.1 M NaOH
Maximum of absorption	346 nm 280 nm		350 nm 278 nm	332 nm 285 nm
E ^{1%} 1cm	268 367		288 344	202 409
3	6470 8860		6950 8300	4860 9870

^{* 1} M HCl + methanol (1+9)

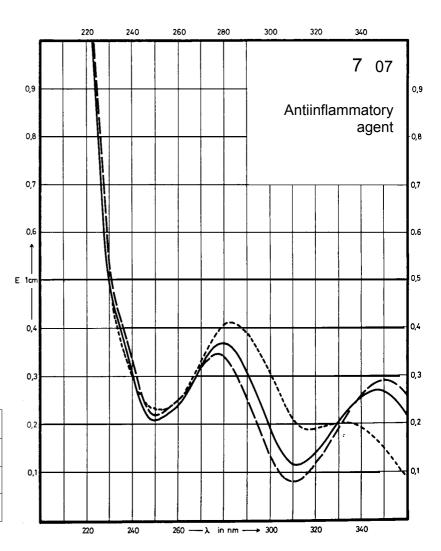

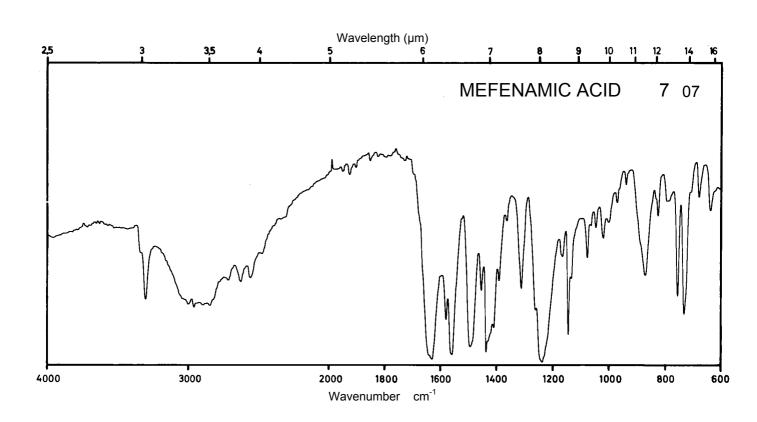


Name METHYLSALICYLATE


COOCH3

 M_{r}


152.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	305 nm 237 nm		303 nm 237 nm	332 nm 243 nm
E ^{1%} 1cm	278 586		245 605	303 416
3	4230 8910		3730 9200	4610 6330

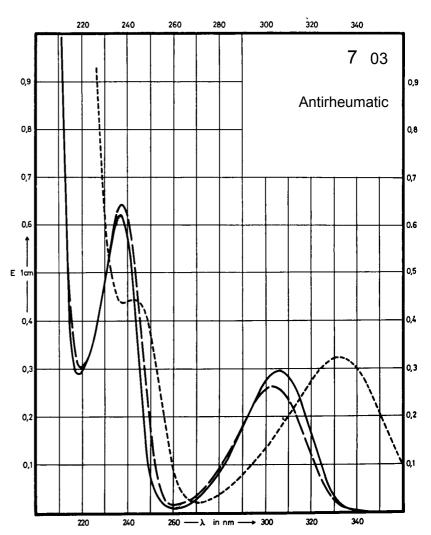

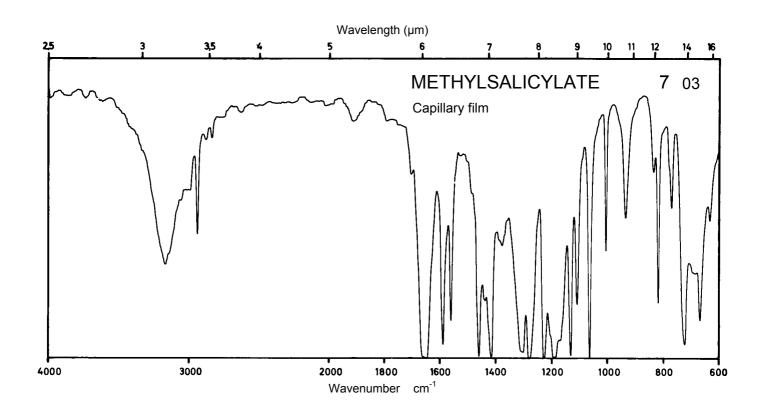




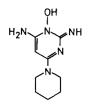


 M_r 209.3


Concentration 0.5 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 262 nm 231 nm	288 nm 263 nm	281 nm 230 nm	288 nm 262 nm
E ^{1%} 1cm	605 544 1650	590 508	547 529	582 512
ε	12700 11400 34500	12300 10600	11400 11100	12200 10700

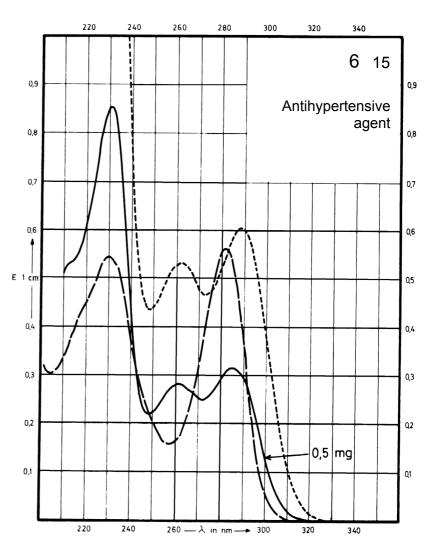

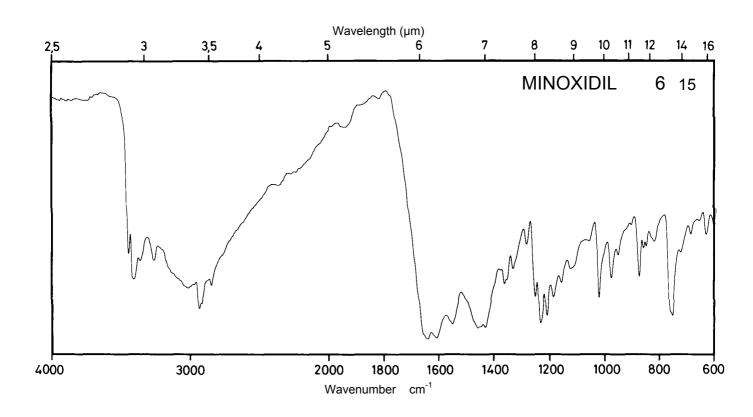



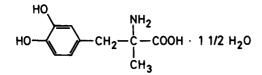



 M_r

Concentration 2 mg / 100 ml


238.2

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of absorption	283 nm		280 nm	
E ^{1%} 1cm	132		119	
3	3150		2830	

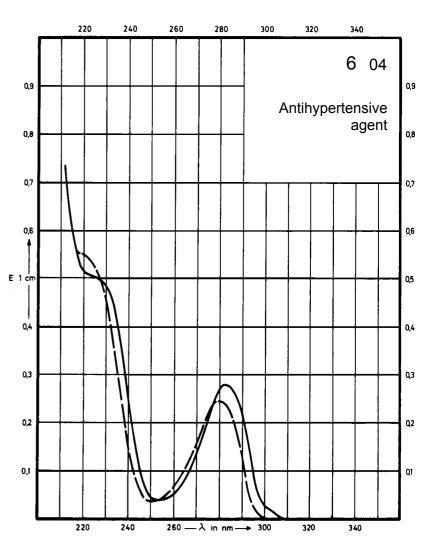

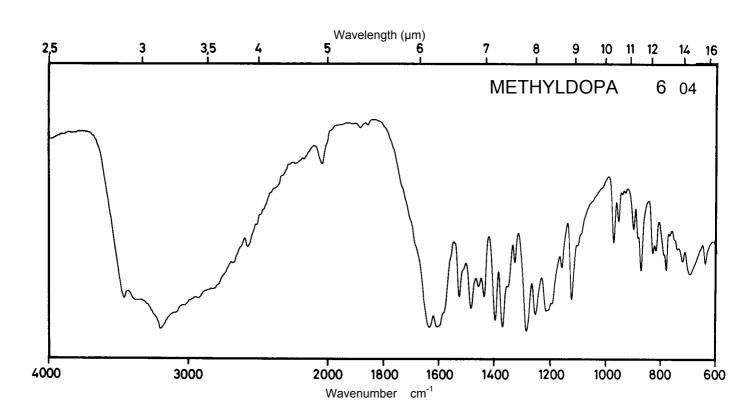




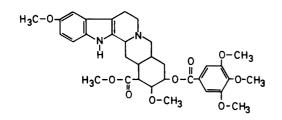


608.7

 M_r

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm		270 nm	284 nm
E ^{1%} 1cm	297		307	163
ε	18080		18690	9920

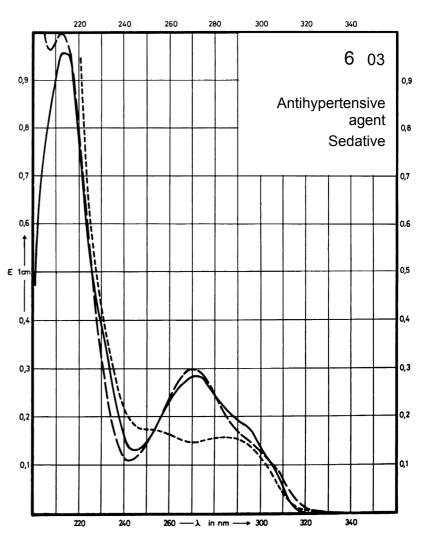


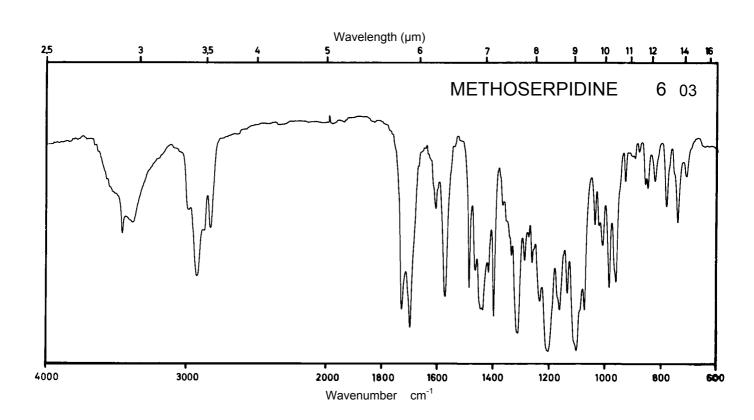



Name METFORMIN HYDROCHLORIDE


 M_r 165.6


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	236 nm	232 nm		232 nm
E ^{1%} 1cm	908	800		787
3	15030	13250		13030

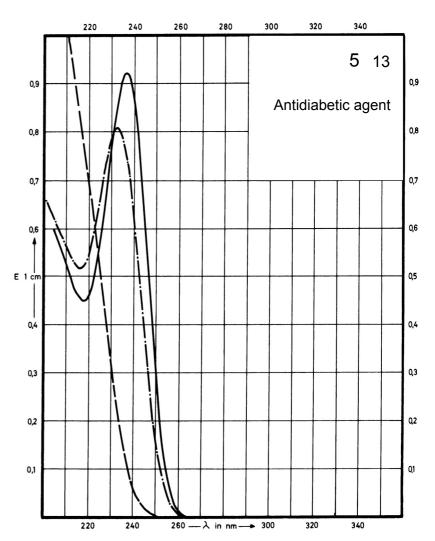

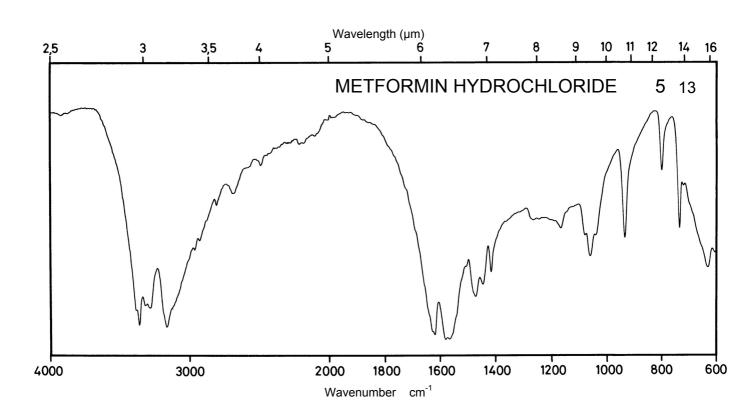

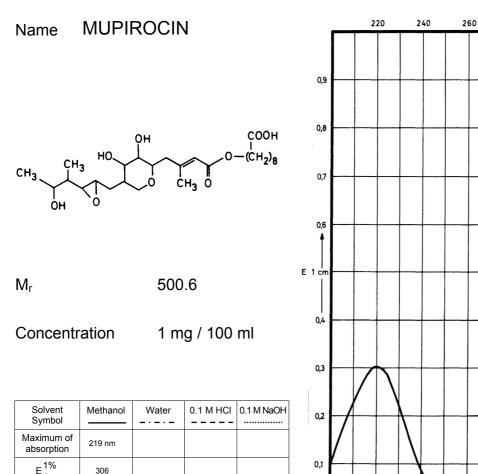

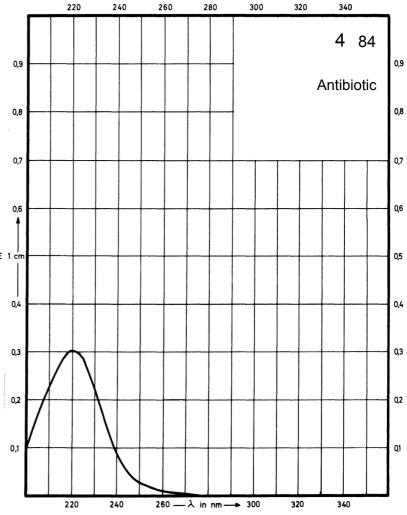


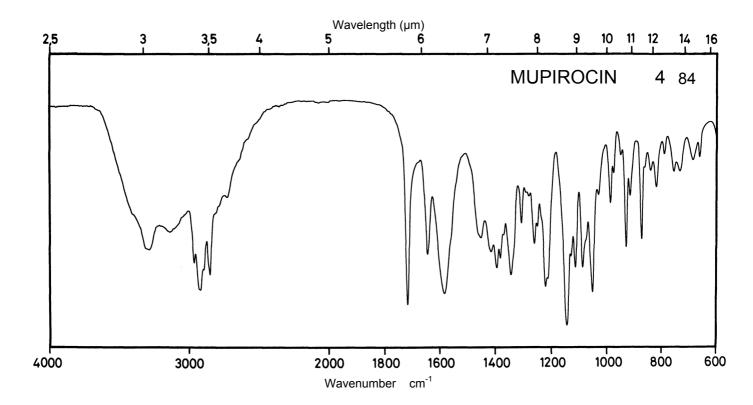

ε

15300

Name MECLOCYCLINE SULFOSALICYLATE

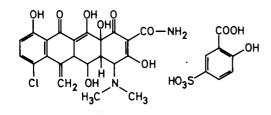



695.1

Concentration


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	348 nm 237 nm		344 nm 238 nm	386 nm 259 nm
E ^{1%} 1cm	180 394		181 398	226 419
ε	12500 27400		12600 27600	15700 29100

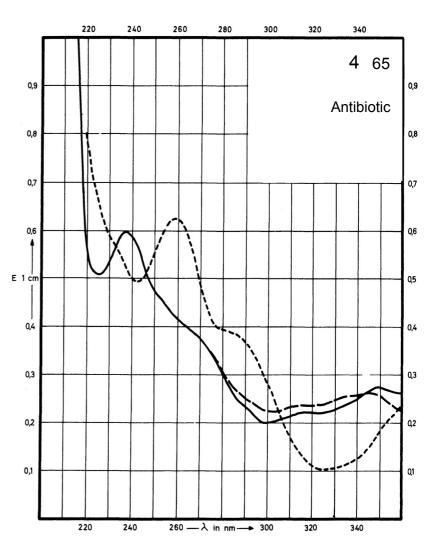

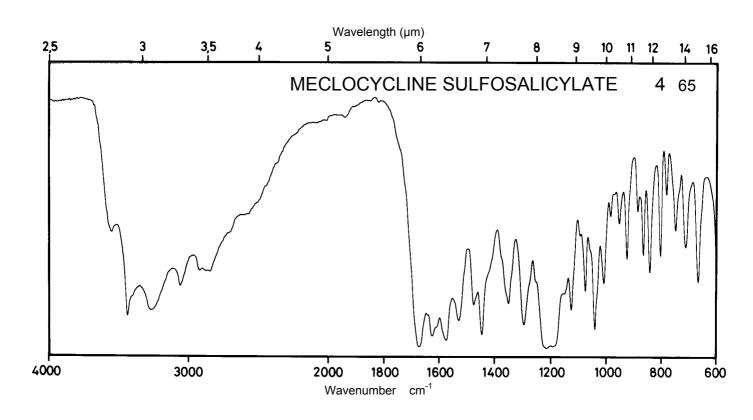


Name MINOCYCLINE HYDROCHLORIDE

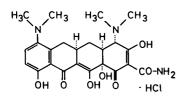


 M_{r}


493.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	348 nm 248 nm		354 nm 265 nm	385 nm 244 nm
E 1% 1cm	256 385		290 370	335 383
ε	12630 19000		14300 18280	16560 18900

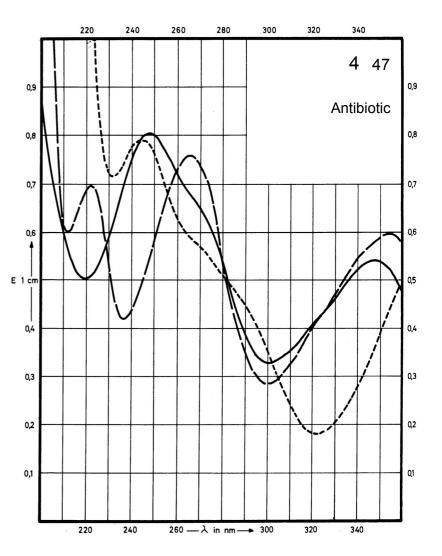

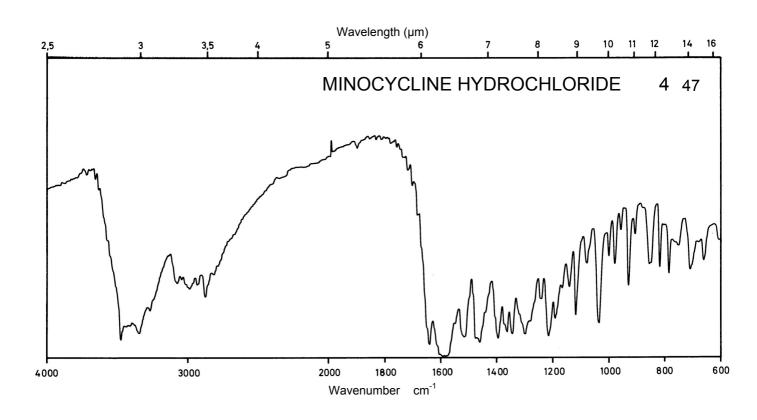




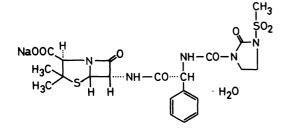


 M_r


579.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 259 nm		Decom- position observed	Decom- position observed
E ^{1%} 1cm	3.17 4.72			
ε	184 274			

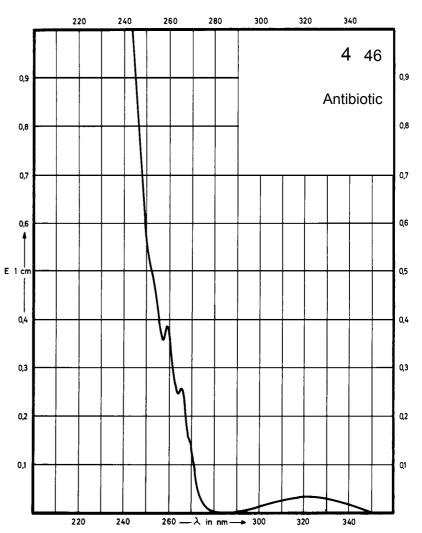

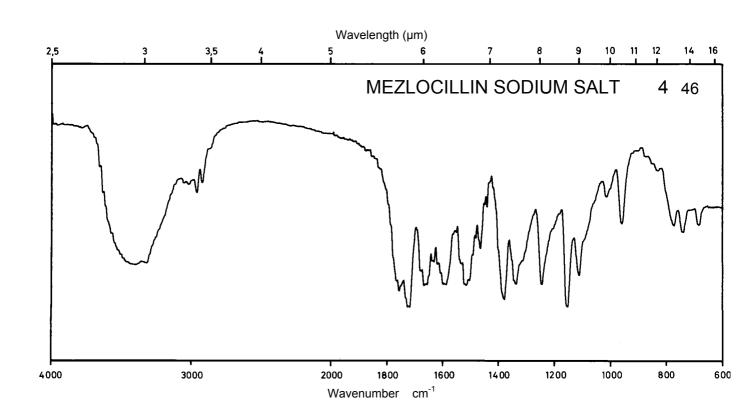


Name METACYCLINE HYDROCHLORIDE

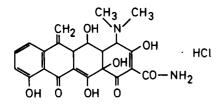


 M_r 478.9

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	349 nm 245 nm	345 nm 282 nm 242 nm	344 nm 245 nm	384 nm 265 nm
E ^{1%} 1cm	314 451	321 328 453	304 457	328 386
ε	15040 21600	15380 15700 21680	14540 21870	15700 18500

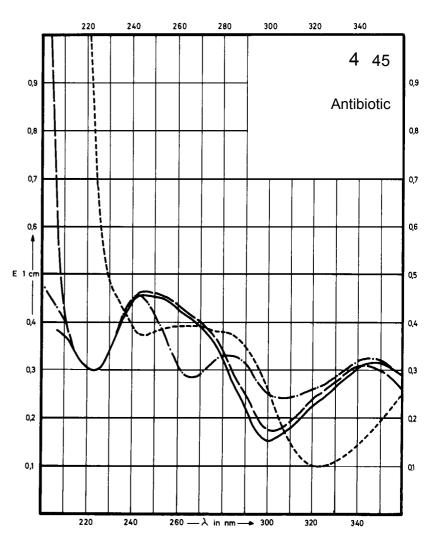

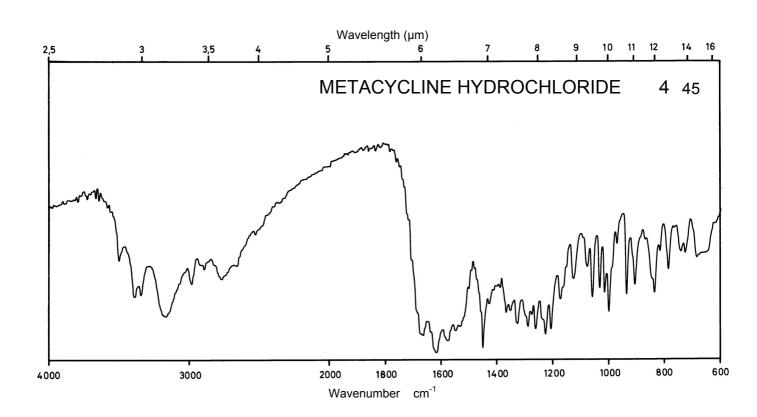


Name METICILLIN SODIUM SALT


 M_{r}

420.4


Concentration 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm	280 nm		
E ^{1%} 1cm	53	55		
ε	2230	2310		

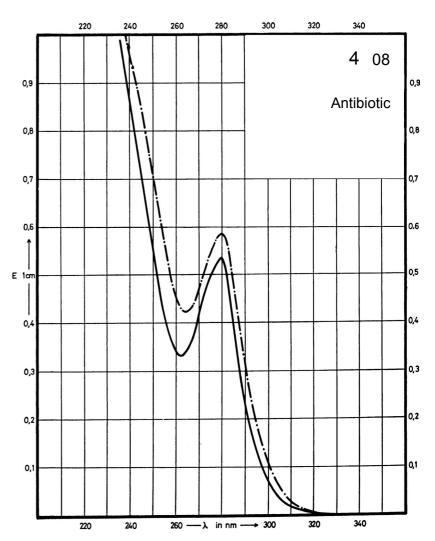

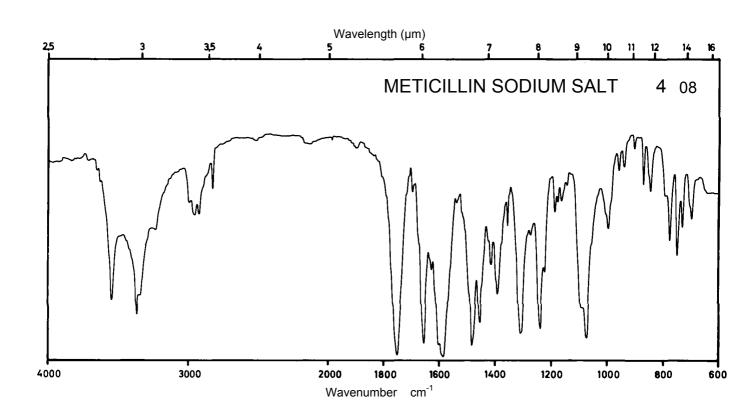


Name MEQUITAZINE

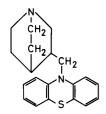


 M_r 322.5


Concentration 0.5 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	305 nm	305 nm	305 nm	306 nm
	253 nm	253 nm	253 nm	254 nm
E 1%	132	119	124	119
1cm	990	928	972	864
ε	4260	3800	4000	3840
	32000	29900	31300	27900

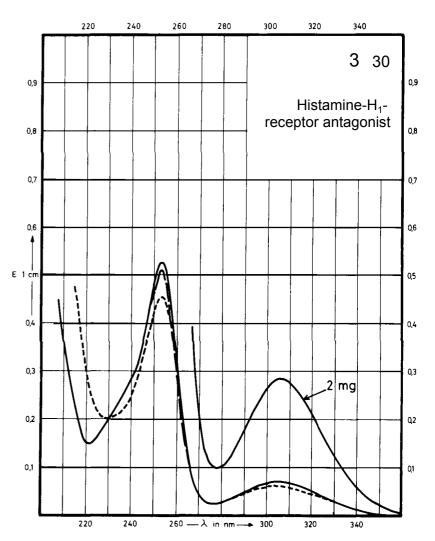

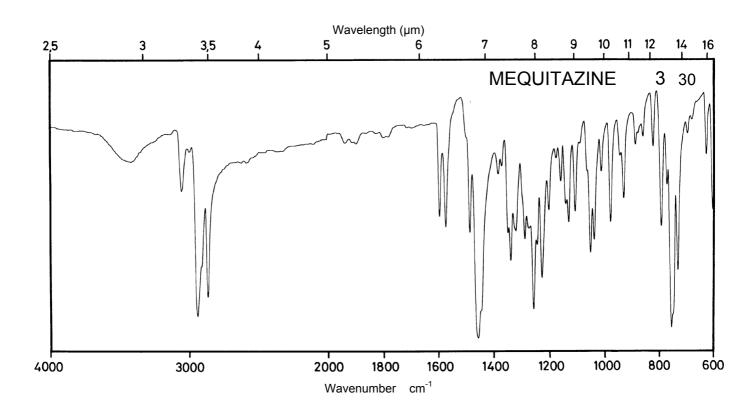




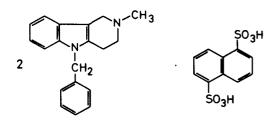


 M_r 841.1


Concentration 2 mg / 100 ml 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	320 nm 315 nm 286 nm		286 nm	287 nm
E ^{1%} 1cm	13 11 273		269	274
ε	1080 930 22980		22640	23020

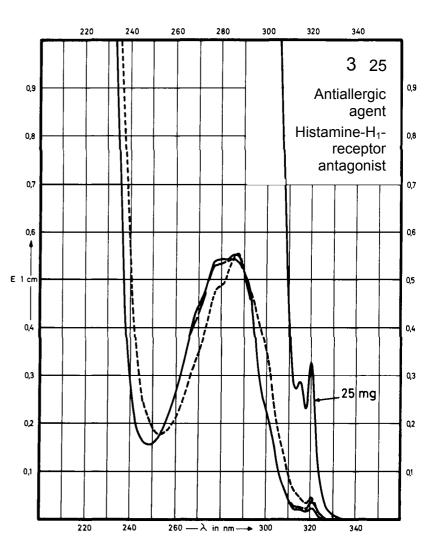

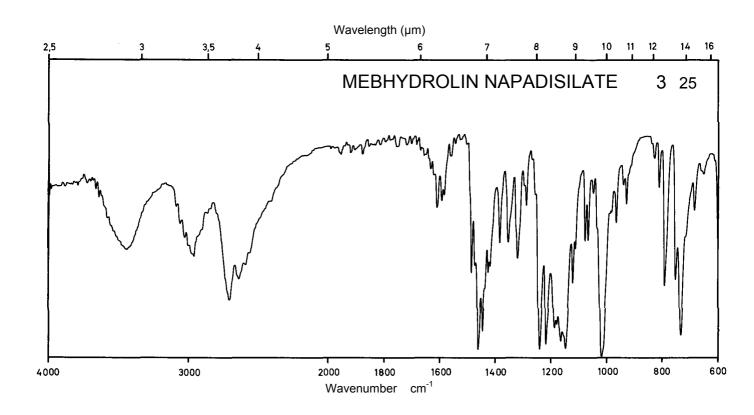



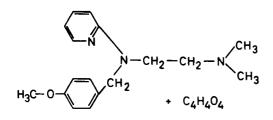



 M_r

401.5


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	308 nm 245 nm		314 nm 236 nm	313 nm 249 nm
E ^{1%} 1cm	122 463		206 440	117 443
ε	4900 18590		8270 17670	4700 17790

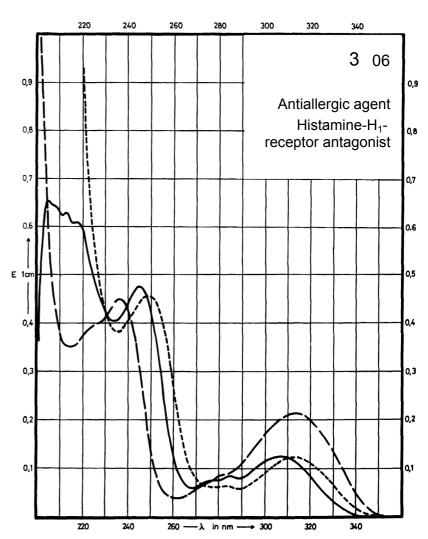

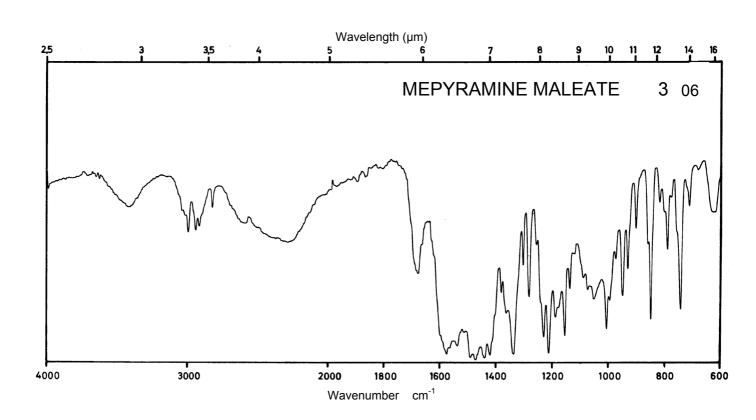




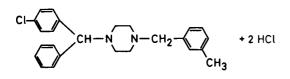


463.9


Concentration 2 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm 229 nm		264 nm 230 nm	
E ^{1%} 1cm	21 309		27 349	
ε	970 14330		1250 16190	

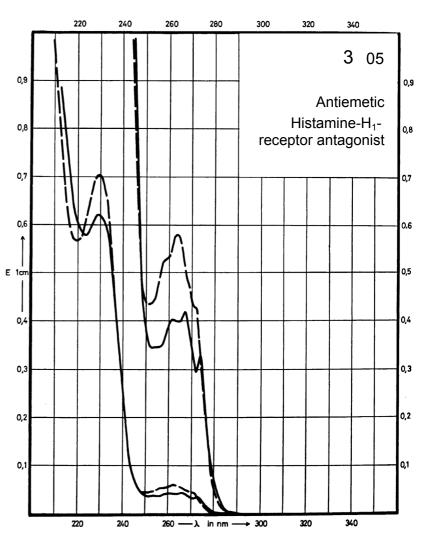

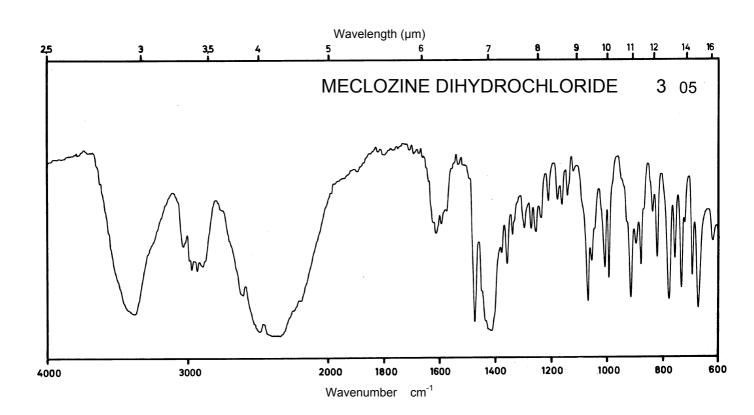




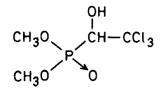


257.4


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
ε				

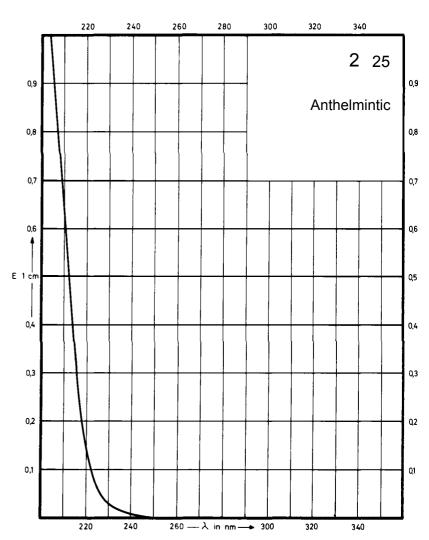

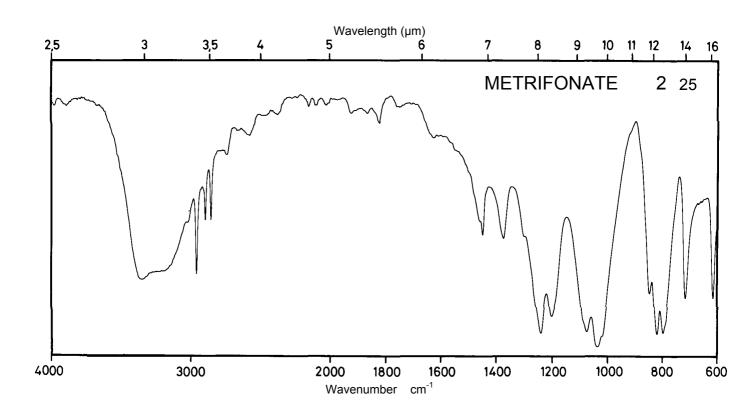

Name MEBENDAZOLE

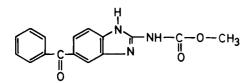


 M_{r}


295.3

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm		286 nm	351 nm 262 nm
E ^{1%} 1cm	509		588	563 812
ε	15030		17360	16630 23980

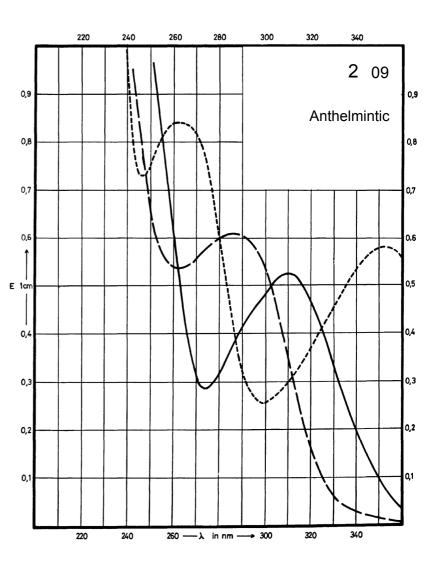

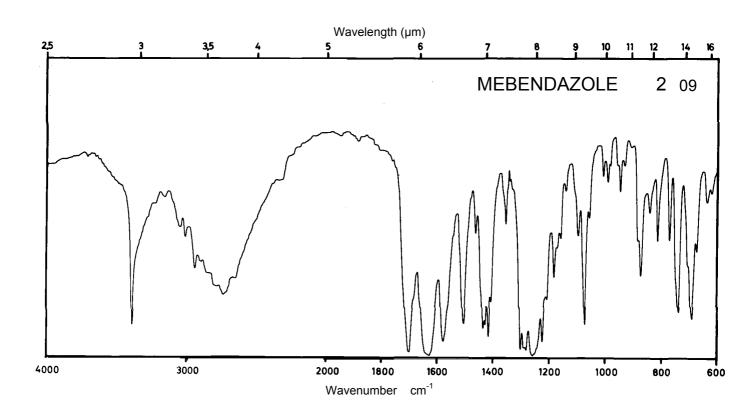

Name MEPTAZINOL HYDROCHLORIDE

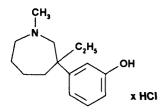



 M_{r}

269.8


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	273 nm 215 nm	273 nm 215 nm	290 nm 237 nm
E ^{1%} 1cm		71 255	70 255	113 323
3		1900 6900	1900 6900	3000 8700

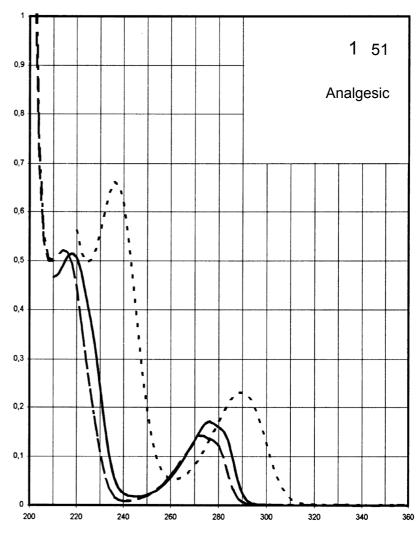

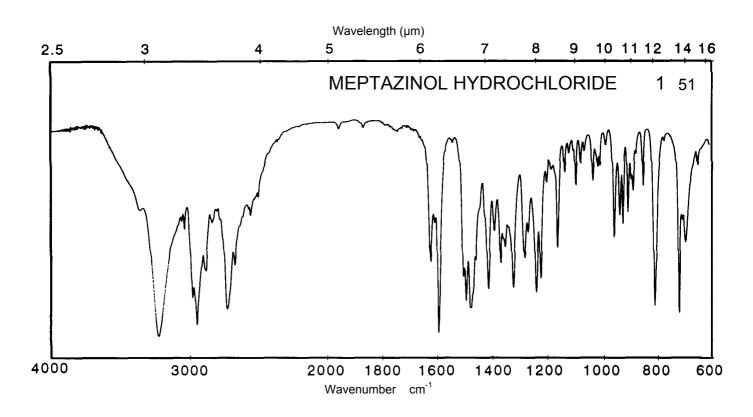




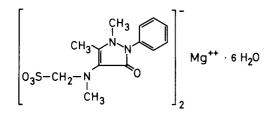


 M_r 753.1


Concentration 3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 234 nm		259 nm	
E ^{1%} _{1cm}	224 265		246	
ε	16830 19960		18530	

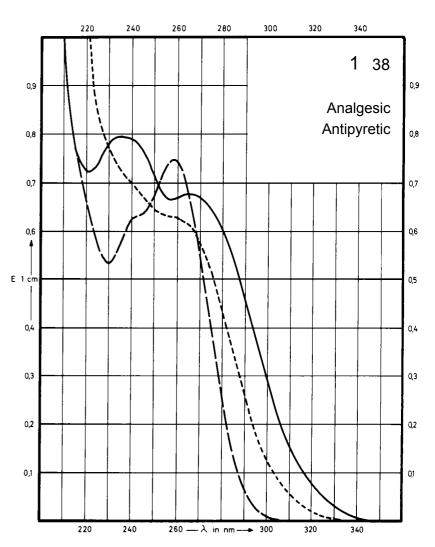

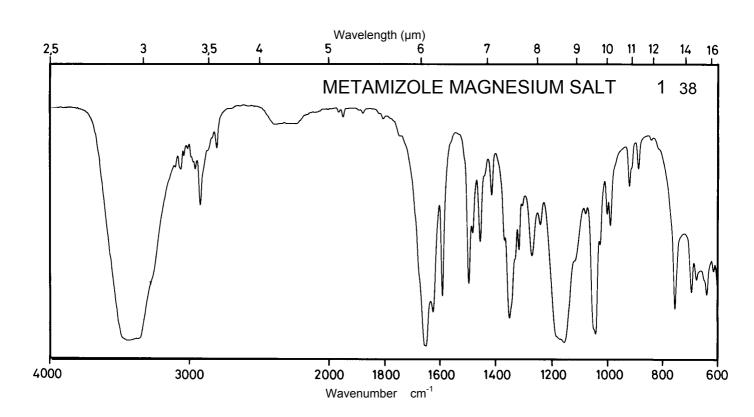




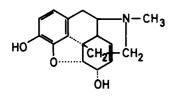


 M_r 285.4

Concentration 14 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm			
E ^{1%} 1cm	53			
3	1500			

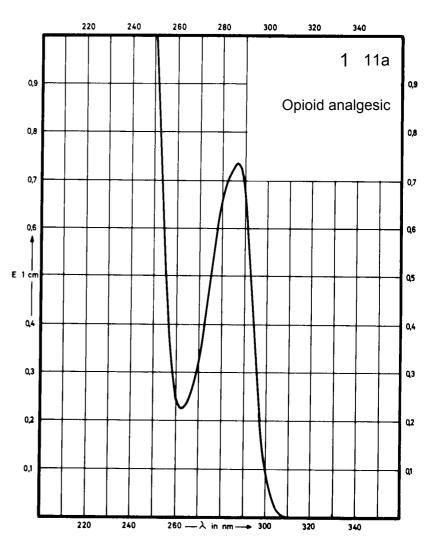

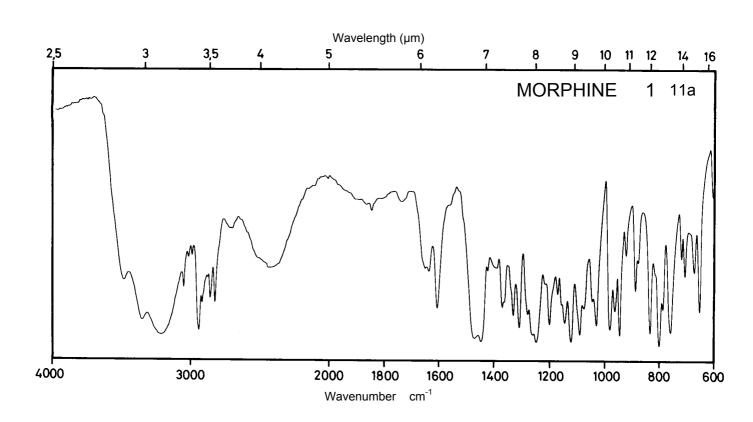


Name MORPHINE HYDROCHLORIDE


 M_{r}

375.9


Concentration 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm		286 nm	297 nm
E ^{1%} 1cm	44		40	68
ε	1650		1500	2560

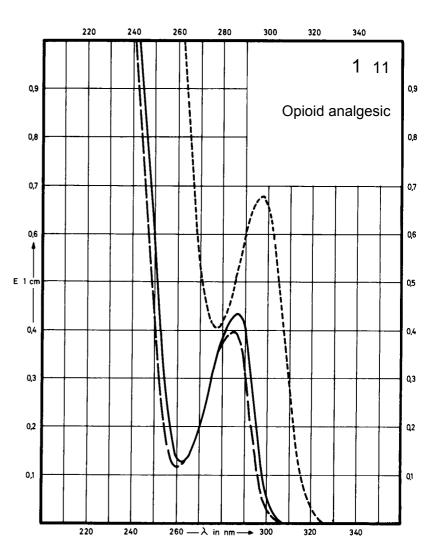

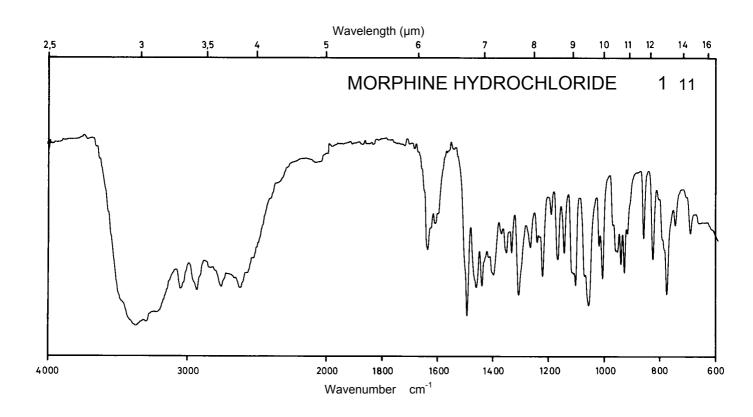




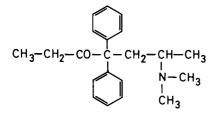


 M_r 309.5

Concentration 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	293 nm			
E ^{1%} 1cm	13.8			
3	428			

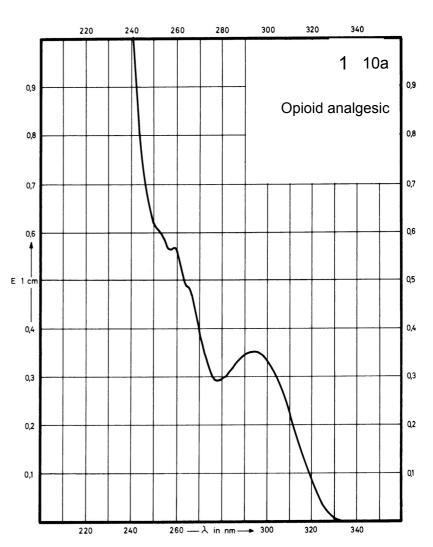

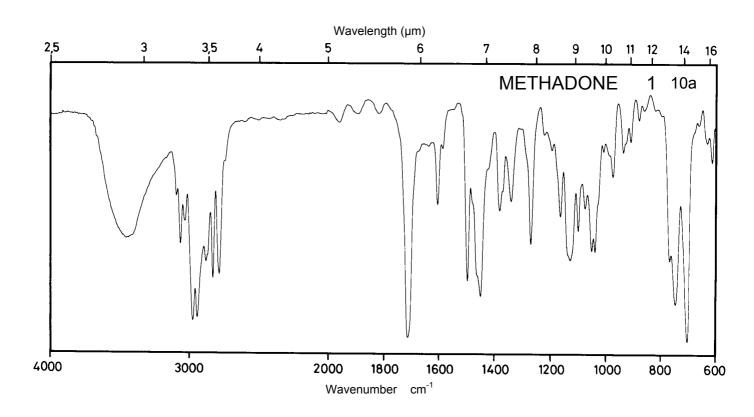


Name METHADONE HYDROCHLORIDE


 M_{r}

345.9


Concentration 50 mg / 100ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	293 nm 258 nm		292 nm 258 nm	
E ^{1%} 1cm	12.3 13.0		16.5 16.4	
3	425 450		570 565	

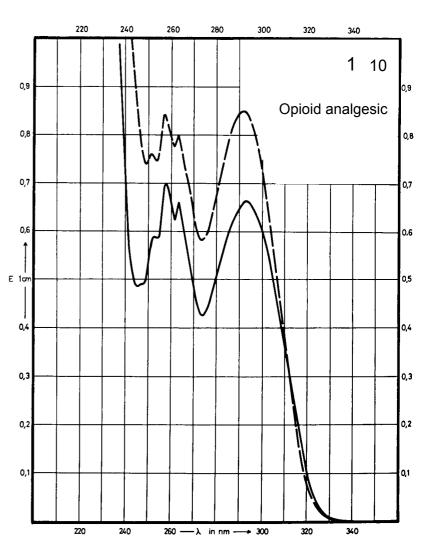

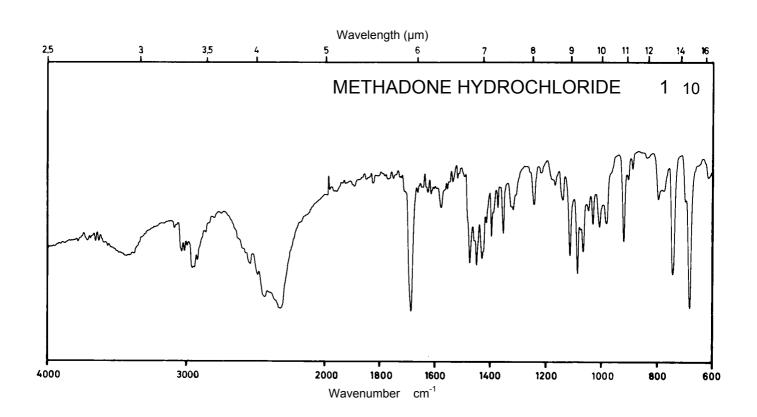

Name METAMIZOLE SODIUM

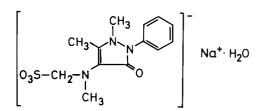



 M_{r}

351.4

Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 234 nm		259 nm	
E ^{1%} _{1cm}	240 288		270	
ε	8440 10120		9490	

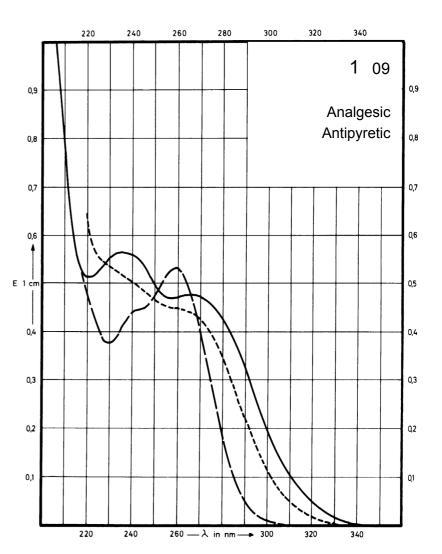

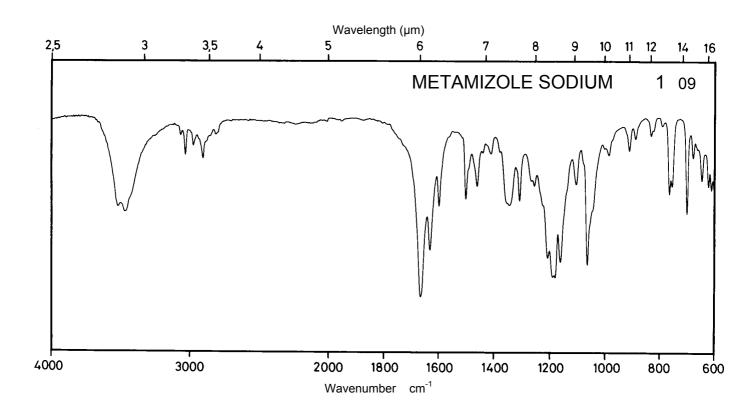


Name METHYLEPHEDRINE HYDROCHLORIDE


 M_{r}

215.7


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm	263 nm 256 nm 251 nm	263 nm 257 nm 251 nm	257 nm 252 nm
E ^{1%} 1cm	6.9 8.7 6.7	7.0 8.9 7.0	7.1 9.1 7.1	9.9 9.0
3	149 188 145	151 192 151	153 196 153	214 194

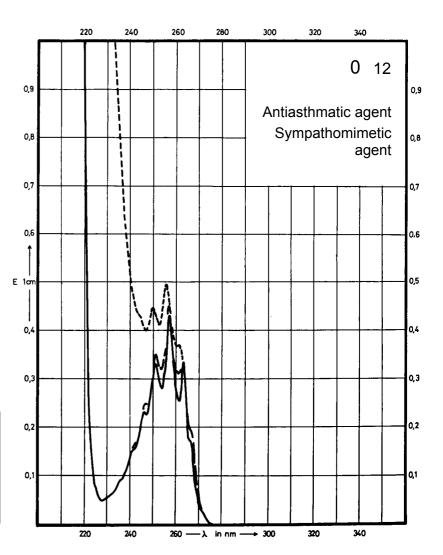

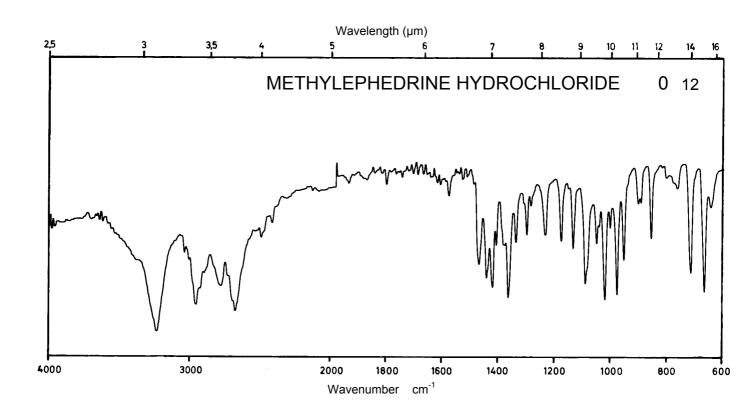




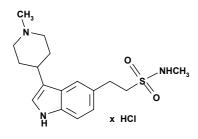


 M_r 371.5


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 225 nm	284 nm 224 nm	284 nm 224 nm	285 nm
E ^{1%} 1cm	163 1088	145 1014	147 1014	147
3	6040 40400	5400 37700	5500 37700	5450

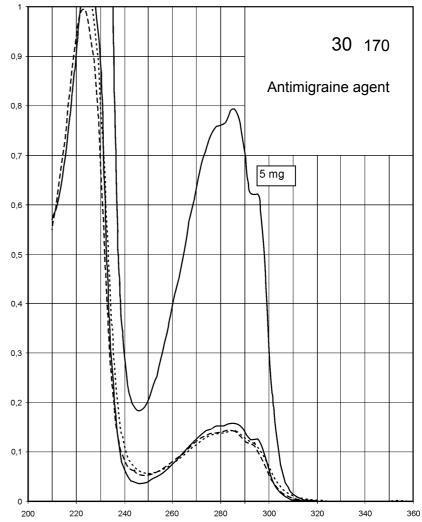

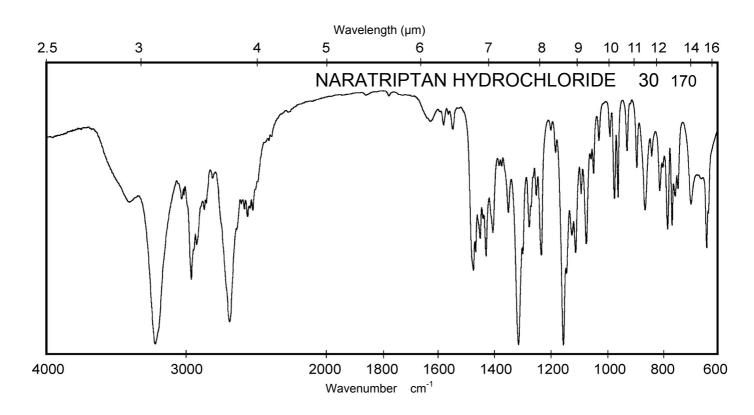



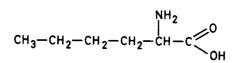



 M_r

131.2


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

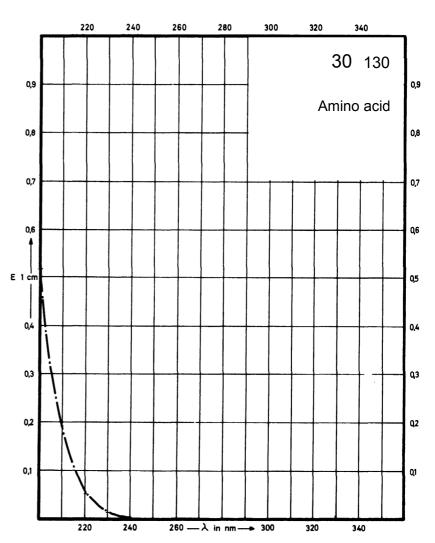

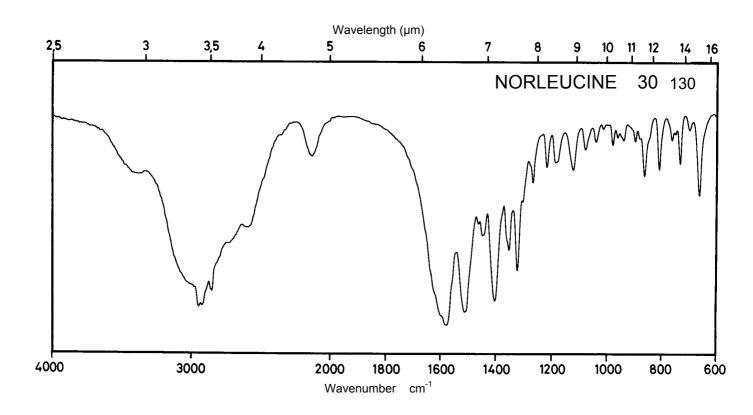




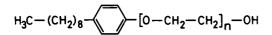


 M_r 616.8


Concentration 29 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 276 nm			
E ^{1%} 1cm	19.8 24.5			
ε	1200 1500			

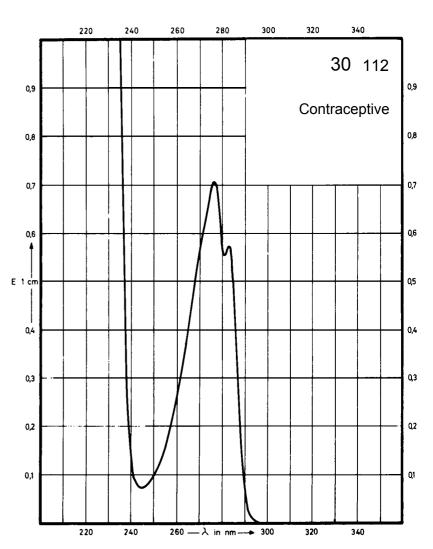

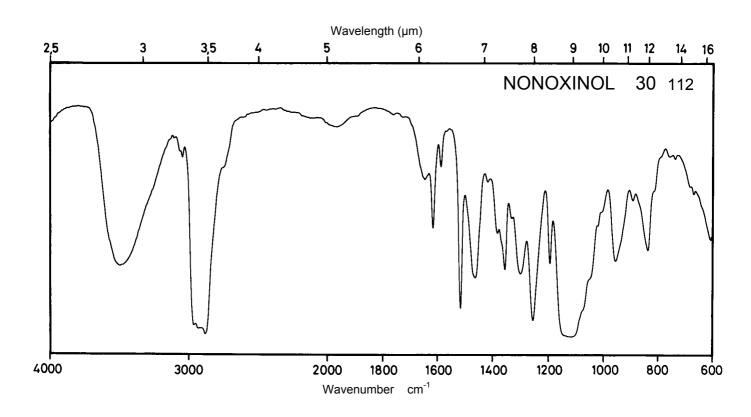




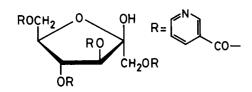


 M_r 600.5


Concentration 1 mg / 100 ml 2 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 269 nm 262 nm 263 nm 261 nm absorption E^{1%} 1cm 201 370 9500 12700 ε 12100 22200

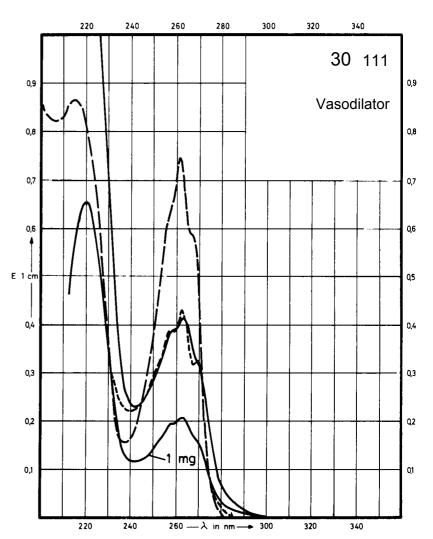

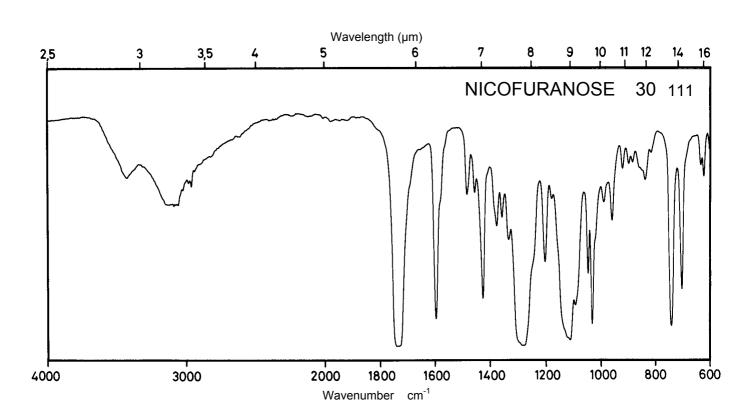

Name NALTREXONE HYDROCHLORIDE

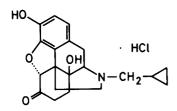



 M_r

377.9


Concentration 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm		281 nm	291 nm
E ^{1%} 1cm	38		33	59
ε	1430		1250	2250

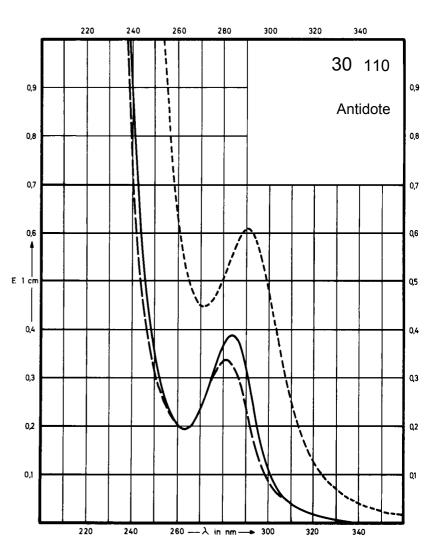

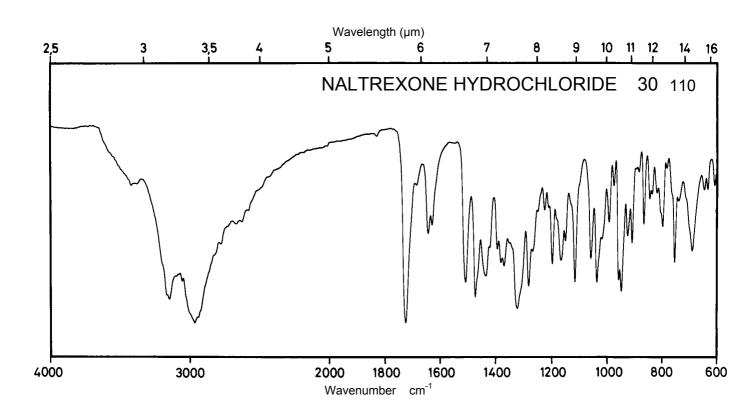




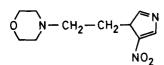


226.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	299 nm		299 nm	305 nm
E ^{1%} 1cm	311		277	353
3	7040		6270	7980

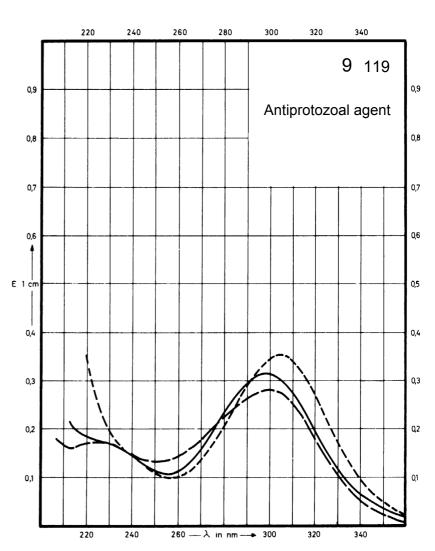

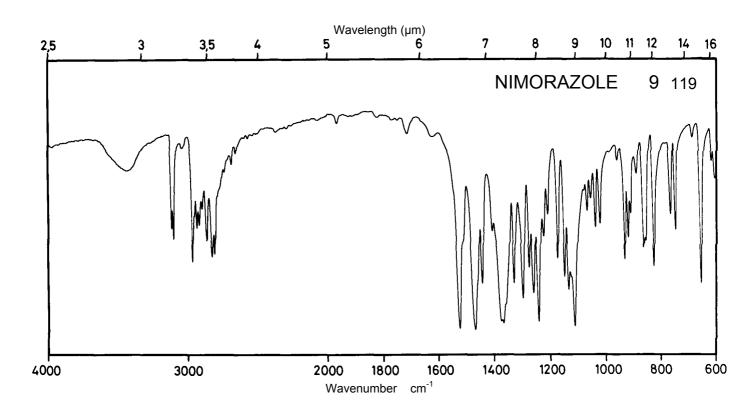




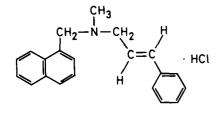


323.9


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	254 nm		254 nm	
E ^{1%} 1cm	675		657	
3	21900		21300	

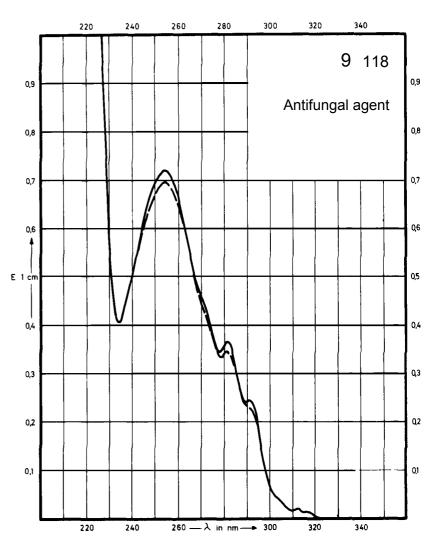

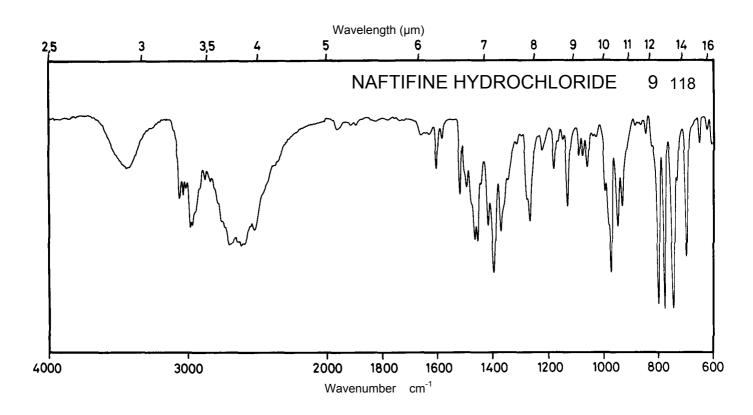




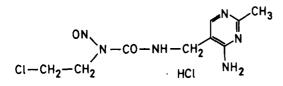


309.2

Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm		245 nm	302 nm
E ^{1%} 1cm	475		498	271
3	14700		15400	8400

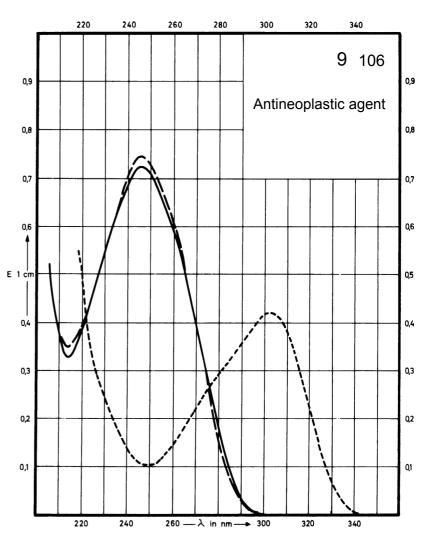

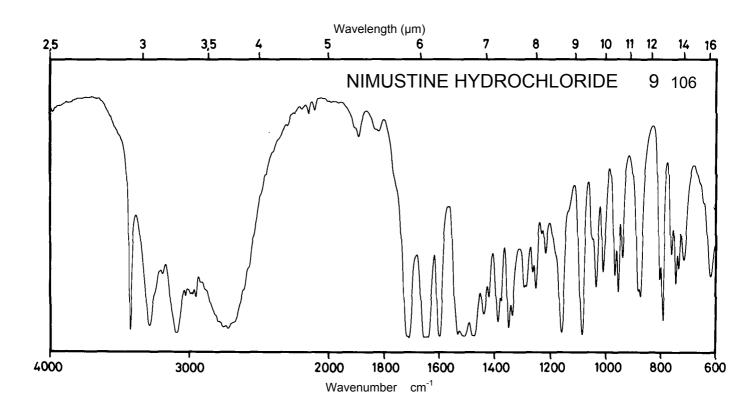

Name NITREFAZOLE


 $\bigvee_{NO_2}^{NO_2} CH_3$


 M_r

248.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	294 nm		305 nm	Decom- position observed
E ^{1%} 1cm	575		570	
ε	14300		14150	

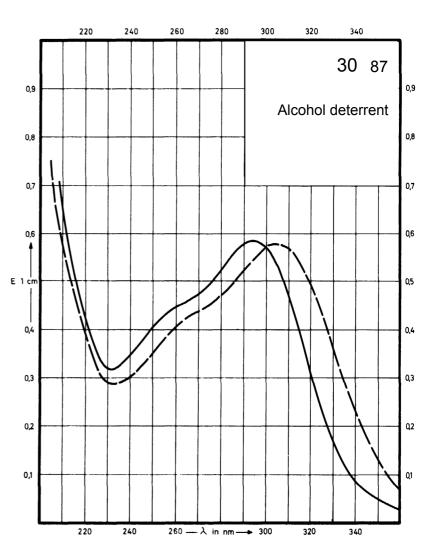

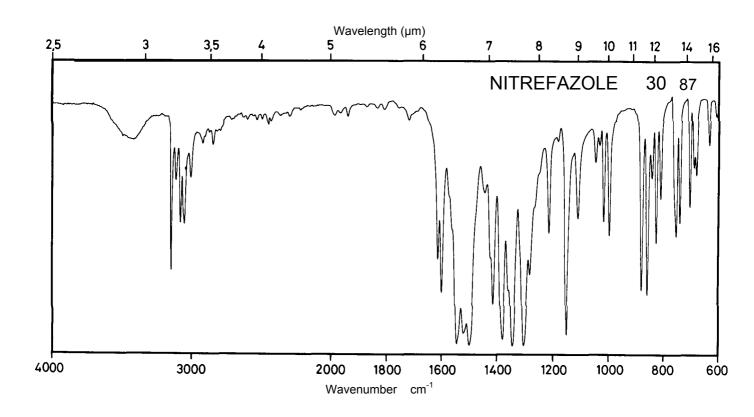


Name NALOXONE HYDROCHLORIDE

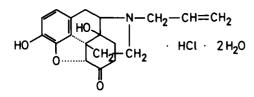


 M_{r} 399.9


Concentration 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm		281 nm	292 nm
E ^{1%} 1cm	35.6		31.4	56.4
ε	1420		1260	2250

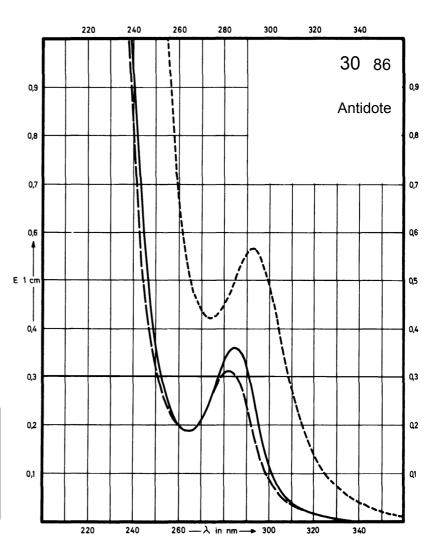

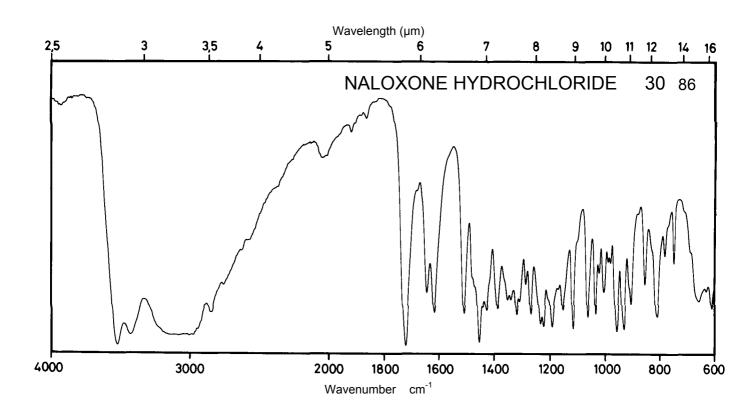


Name NALBUPHINE HYDROCHLORIDE

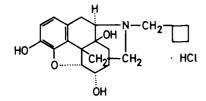


 M_r 393.9

Concentration 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm		284 nm	297 nm
E ^{1%}	46.2		41.1	65.4
ε	1800		1600	2600

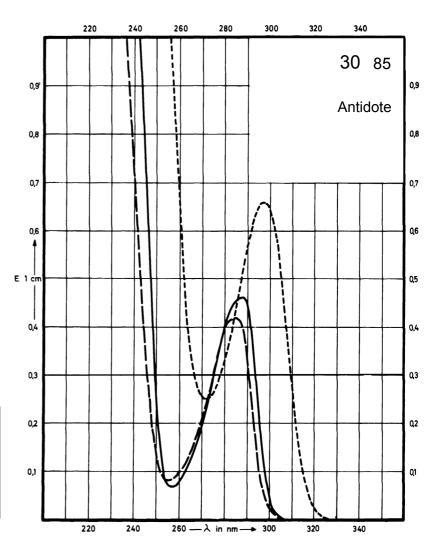

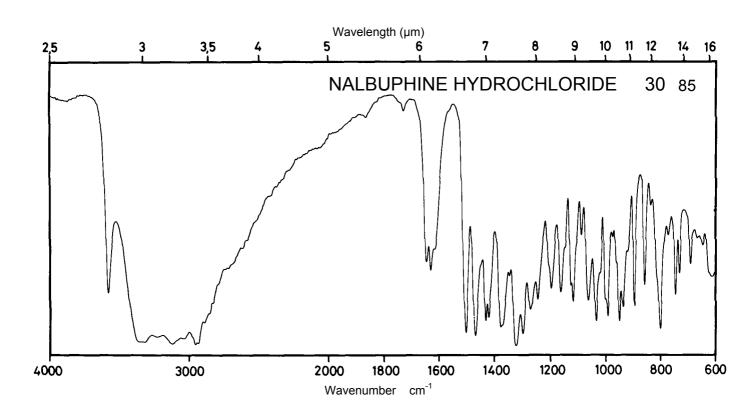


Name NORPSEUDO-EPHEDRINE HYDROCHLORIDE


 M_{r}

187.7

Concentration 70 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm	263 nm 257 nm 252 nm		
E ^{1%} 1cm	7.80 9.82 7.43	8.04 10.2 7.84		
3	146 184 140	150 191 147		

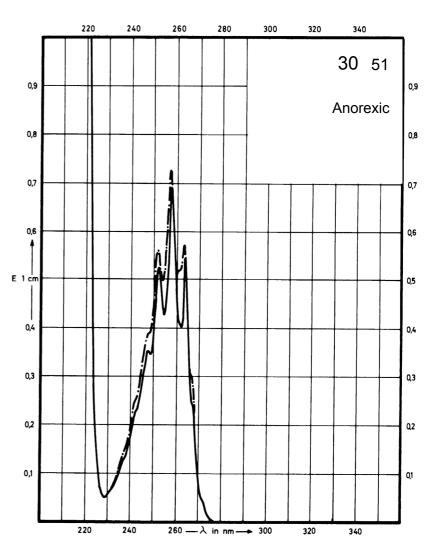

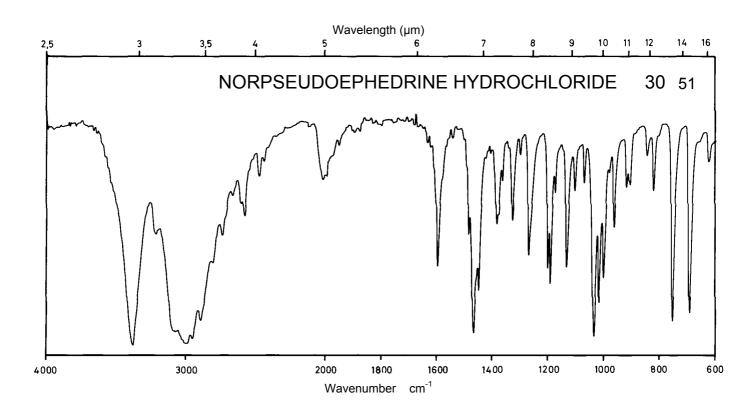


Name NALORPHINE HYDROBROMIDE


HO
$$CH_2$$
 CH_2 CH_2

 M_r 392.3

Concentration 5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm		286 nm	299 nm 251 nm
E ^{1%} 1cm	42		39	66 151
ε	1650		1530	2580 5920

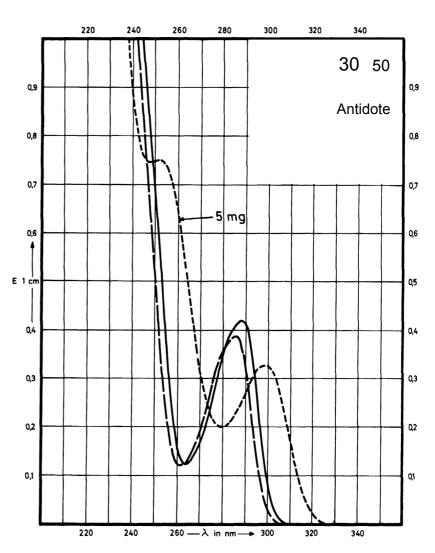

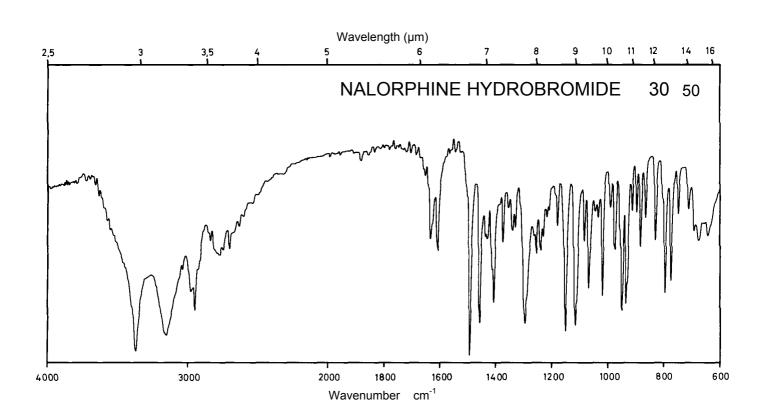







162.2

Concentration


2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm		259 nm	260 nm
E ^{1%} 1cm	181		349	201
ε	2940		5660	3260

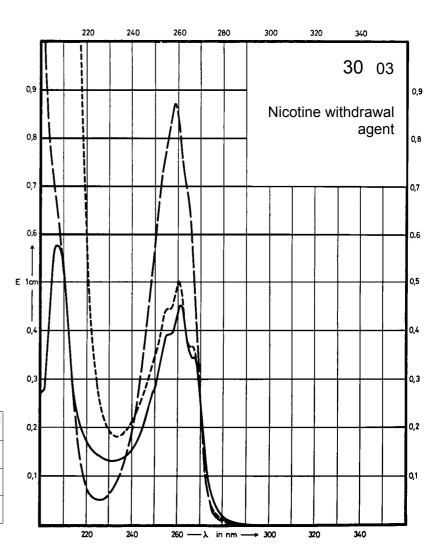

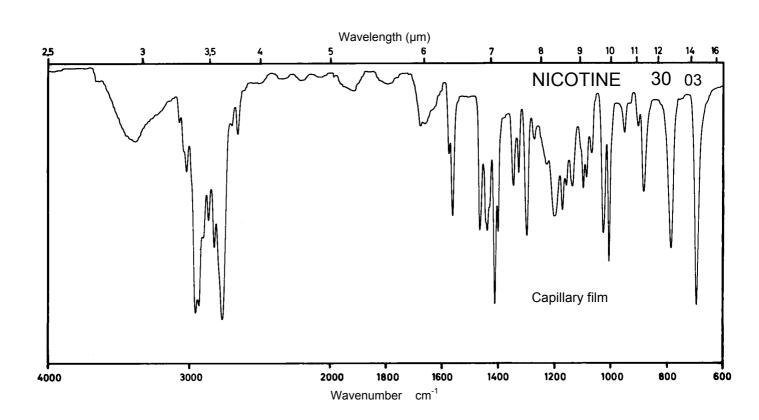




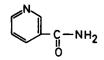


122.1


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm		260 nm	262 nm
E ^{1%} 1cm	221		410	231
ε	2700		5010	2820

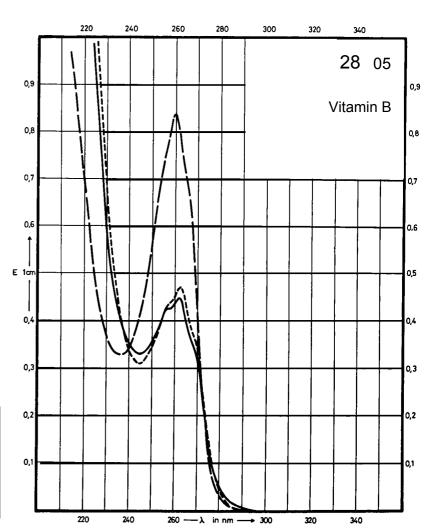

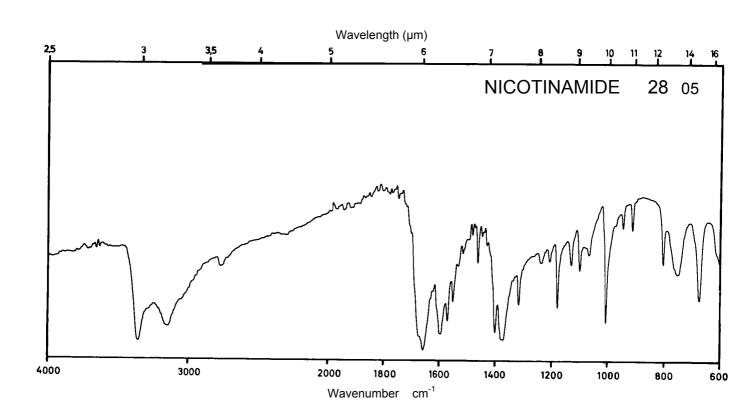




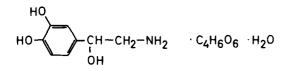


337.3


Concentration 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm		280 nm	Decom- position observed
E ^{1%} 1cm	91		80	
3	3070		2700	

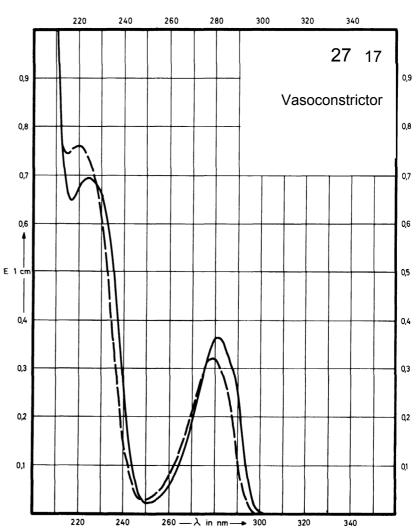

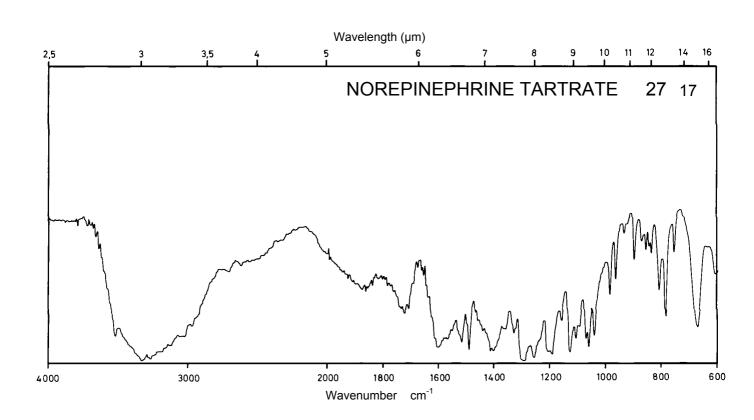




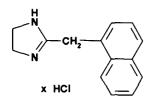


246.7


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 281 nm 271 nm	280 nm 271 nm	280 nm 271 nm	280 nm 271 nm
E ^{1%} 1cm	203 295 247	277 235	278 236	261 216
ε	5000 7300 6100	6800 5800	6900 5800	6400 5300

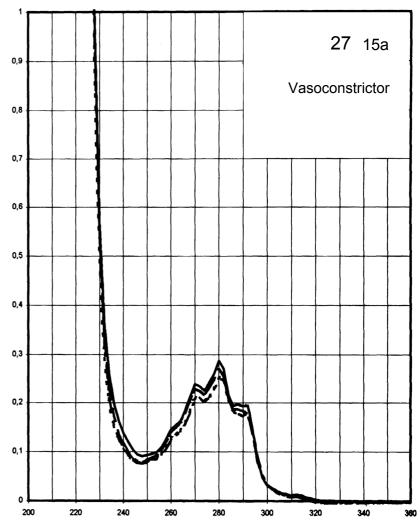

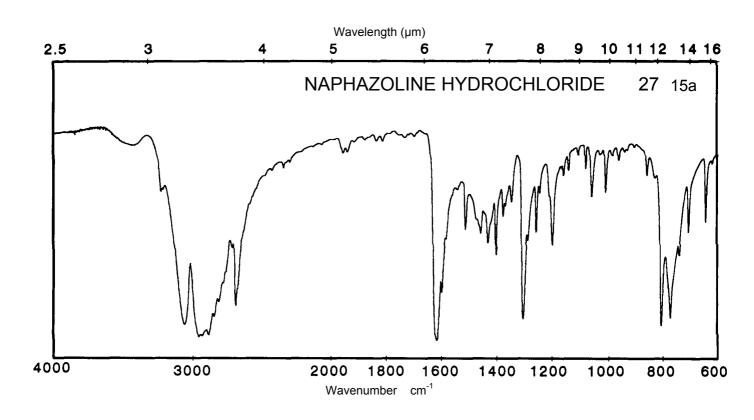

Name NAPHAZOLINE NITRATE

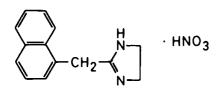



 M_{r}

273.3


Concentration 2 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	313 nm 280 nm 270 nm		291 nm 280 nm 270 nm	281 nm 271 nm
E ^{1%} 1cm	24 265 220		170 253 214	250 211
ε	655 7240 6010		4650 6910 5850	6830 5770

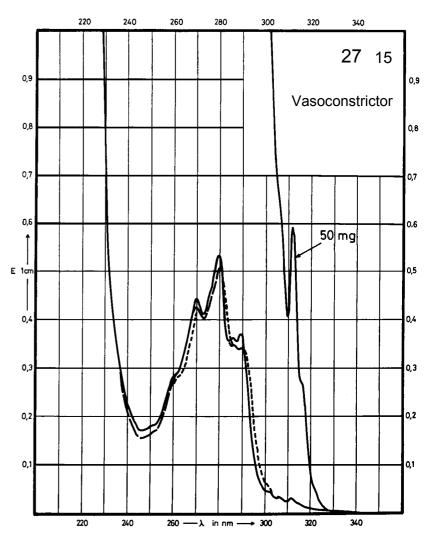

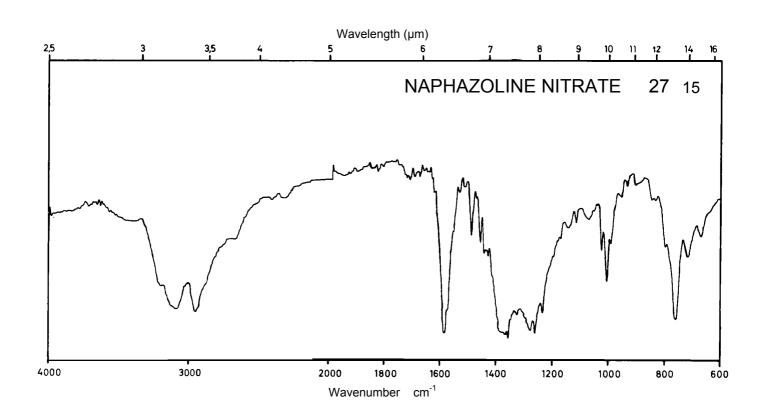


Name NORFENEFRINE HYDROCHLORIDE

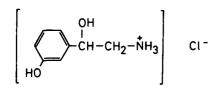


 M_r 189.6


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm		272 nm	290 nm 237 nm
E ^{1%} 1cm	111		96	158 473
ε	2100		1820	3000 8970

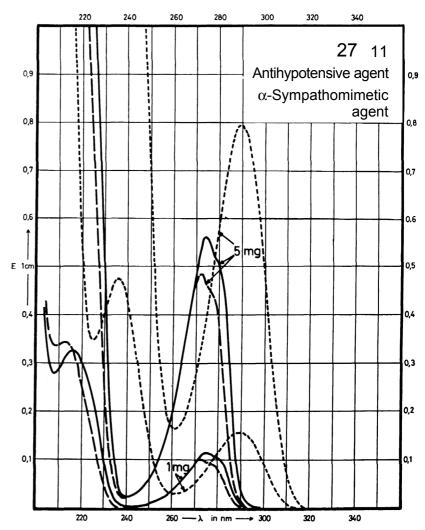

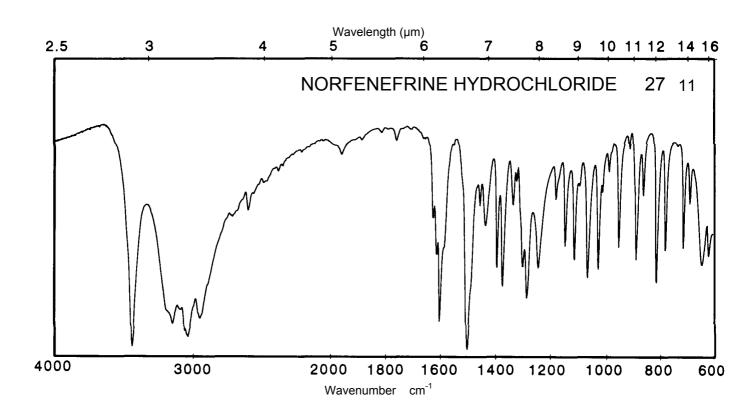




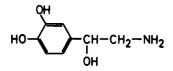


169.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		278 nm	
E ^{1%} 1cm	180		160	
3	3050		2710	

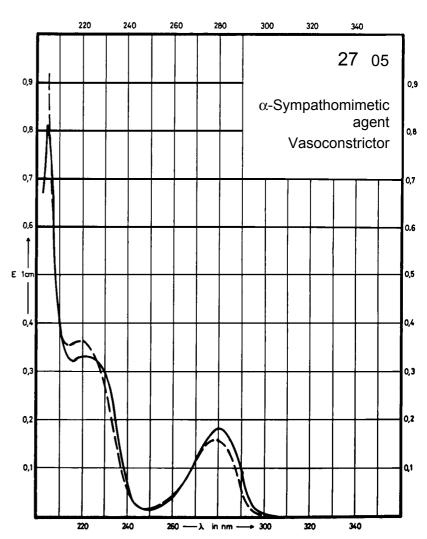

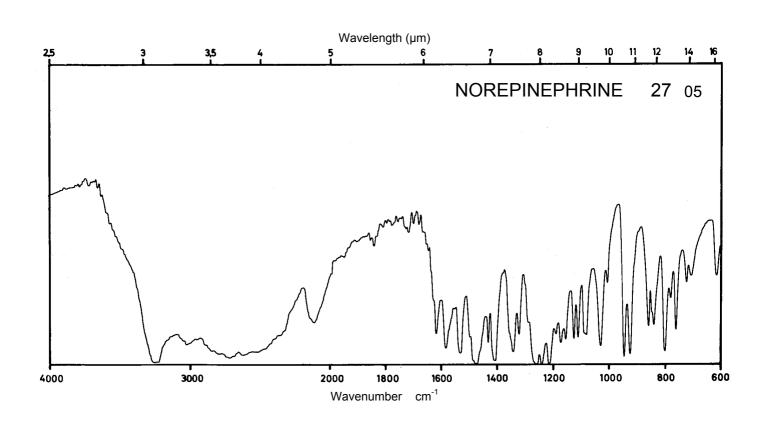




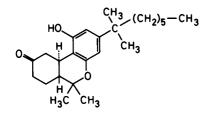


 M_r


372.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm 276 nm			
E ^{1%} 1cm	37 35			
ε	1360 1300			

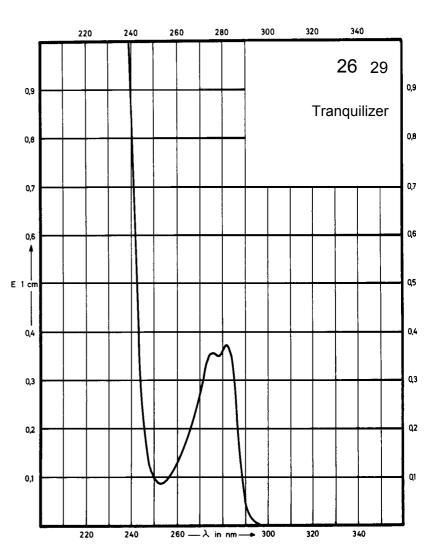

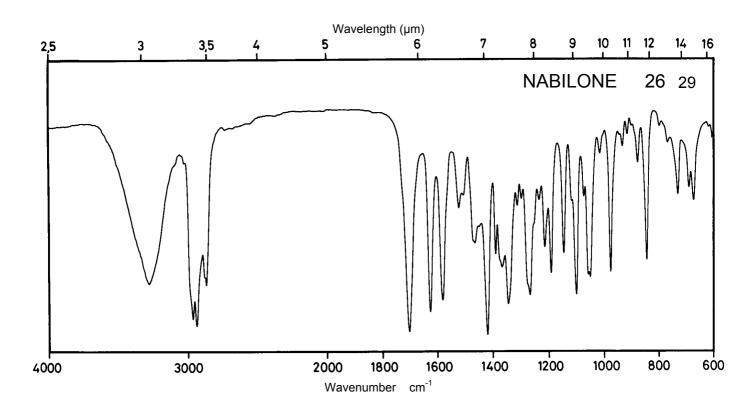





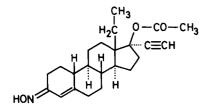


369.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm			
E ^{1%} 1cm	546			
ε	20200			

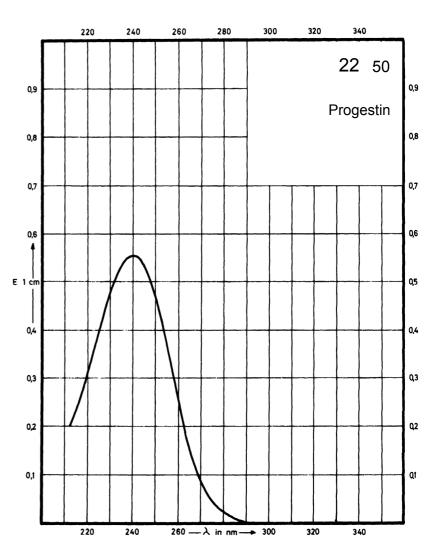

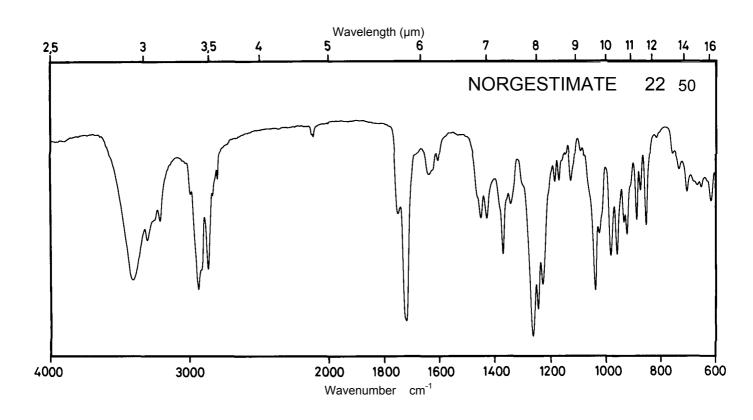





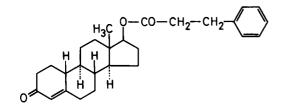


 M_r 406.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm			
E ^{1%} 1cm	432			
ε	17600			

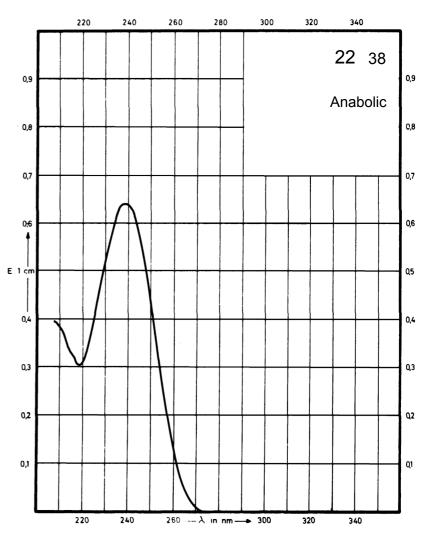

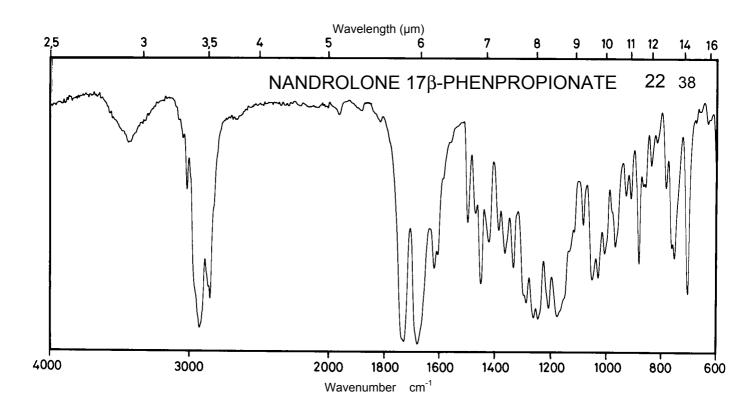



Name NANDROLONE 17β -DECANOATE

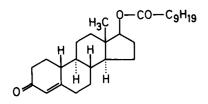


428.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm			
E ^{1%} 1cm	402			
ε	17200			

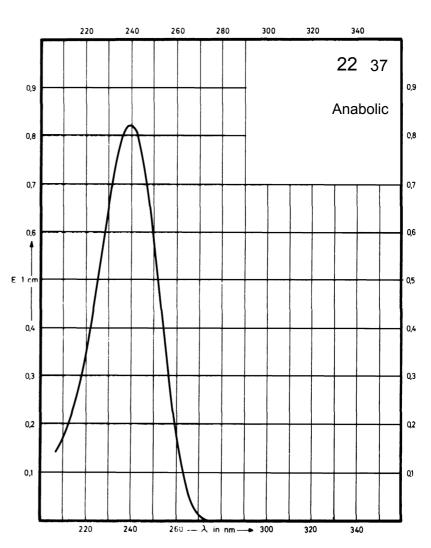

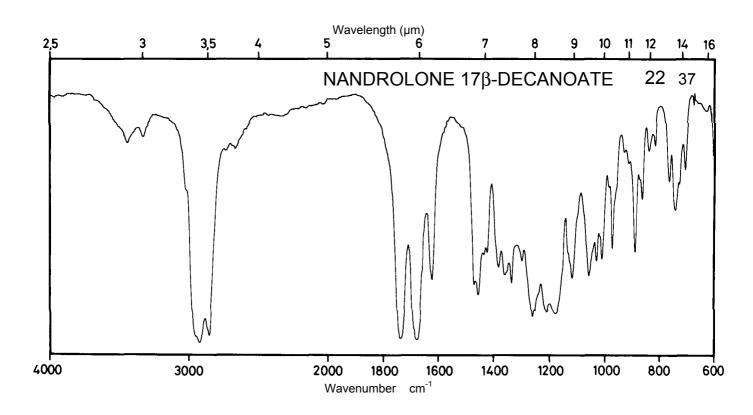





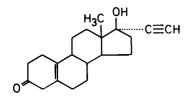


298.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

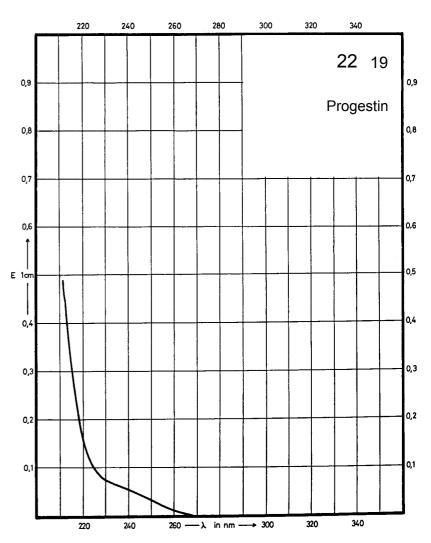

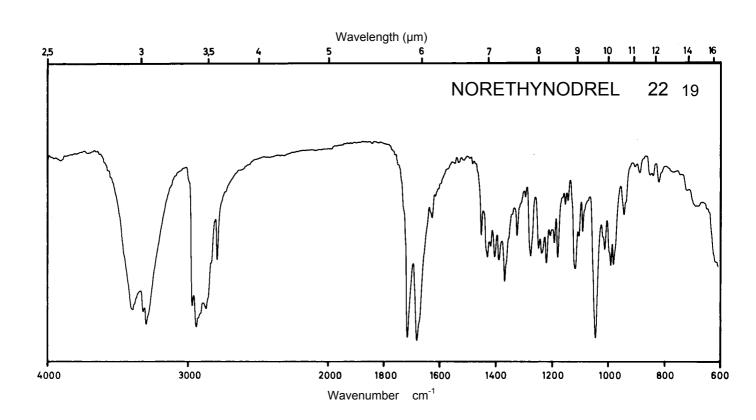


Name NORGESTREL

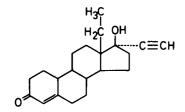


 M_{r}


312.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm		247 nm	247 nm
E ^{1%} 1cm	554		543	527
ε	17310		16970	16480

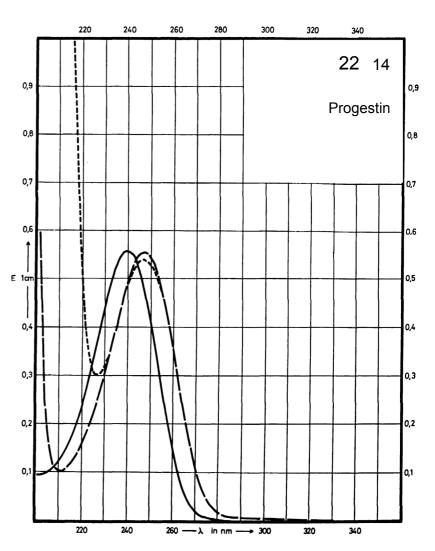

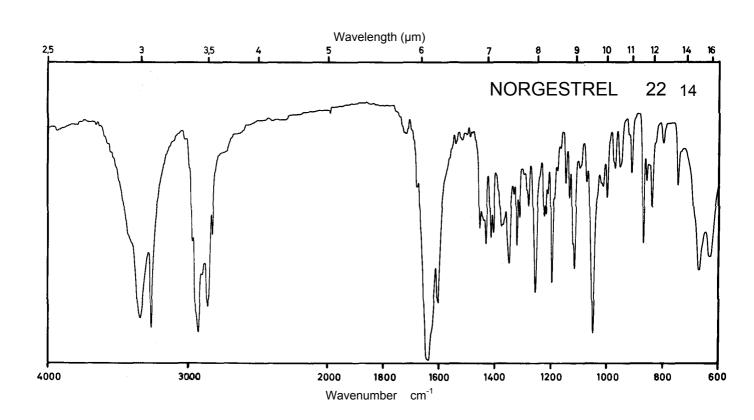


Name NORETHISTERONE

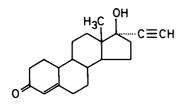


 M_{r}


298.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm		247 nm	247 nm
E ^{1%} 1cm	573		574	564
3	17100		17130	16840

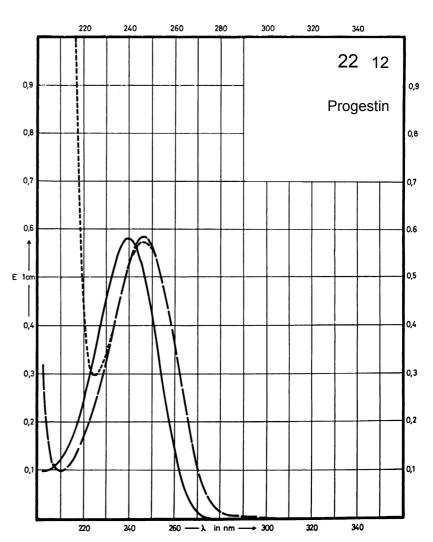

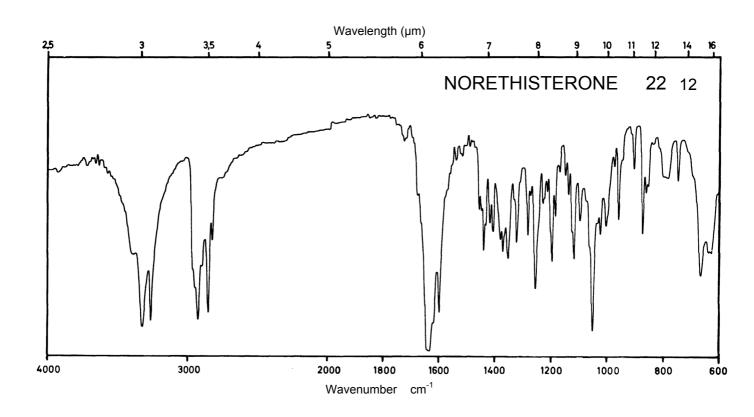


Name NORETHISTERONE 17β-ACETATE

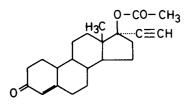


 M_{r}


340.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm		247 nm	247 nm
E ^{1%} 1cm	510		495	495
ε	17390		16900	16900

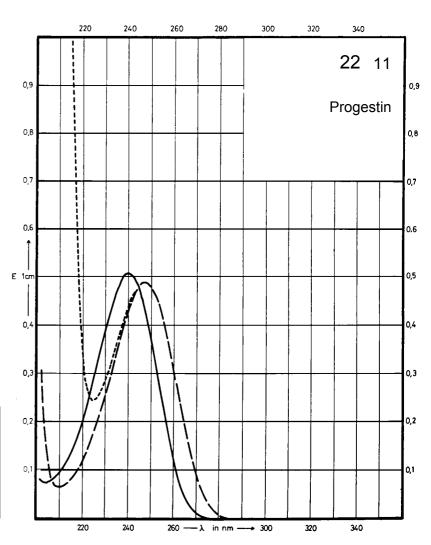

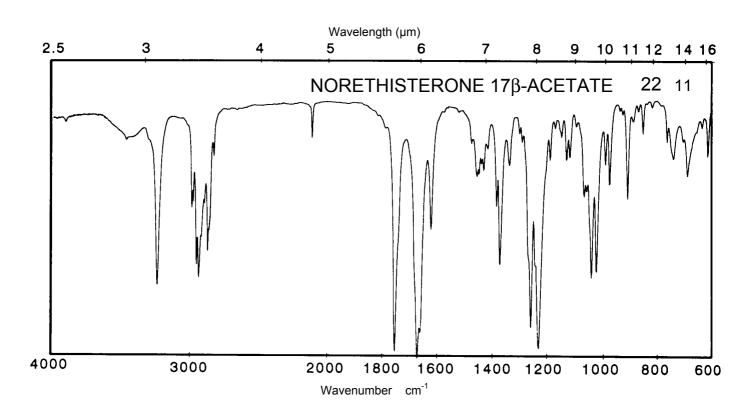

Name NORGESTRIENONE

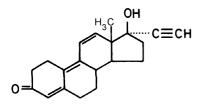


 M_{r}


294.4

Concentration


0.9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	340 nm		350 nm 237 nm	350 nm
E ^{1%} 1cm	987		946 258	938
ε	29060		27850 7600	27610

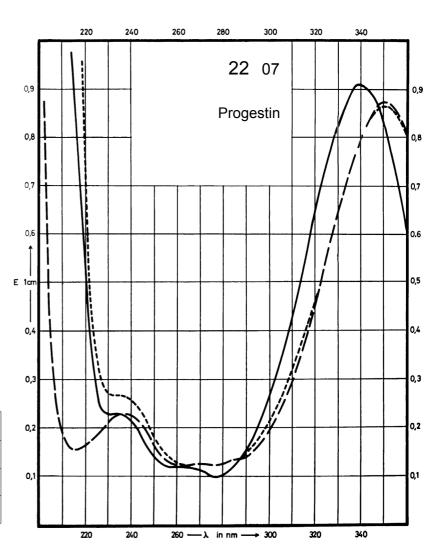

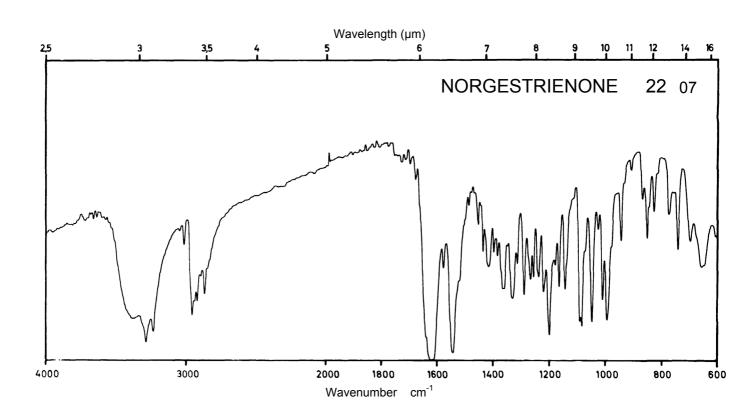


Name NOXIPTILINE HYDROCHLORIDE

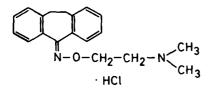


 M_{r}

330.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	250 nm		250 nm	251 nm
E ^{1%} 1cm	410		408	389
ε	13550		13490	12880

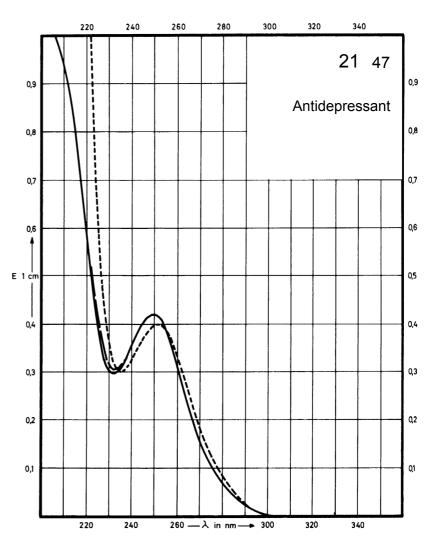

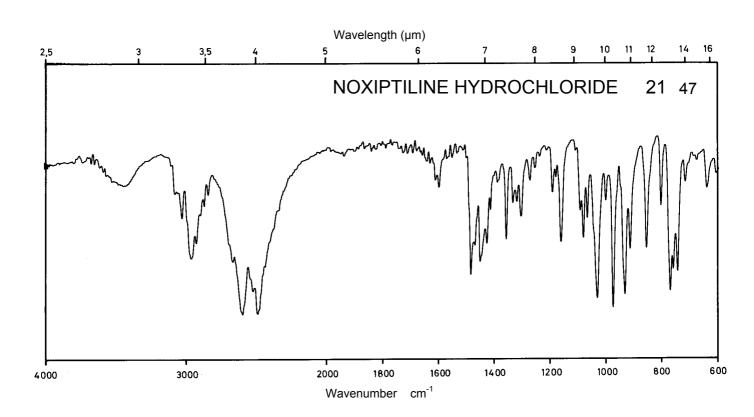


Name NOMIFENSINE MALEATE


 M_{r}


354.4


Solvent Symbol	Methanol	Water	0.1 M HCI	NH ₃ *
Maximum of absorption	293 nm			286 nm
E ^{1%} _{1cm}	78			58
ε	2760			2060

^{* 2.5} per cent NH₃ in methanol

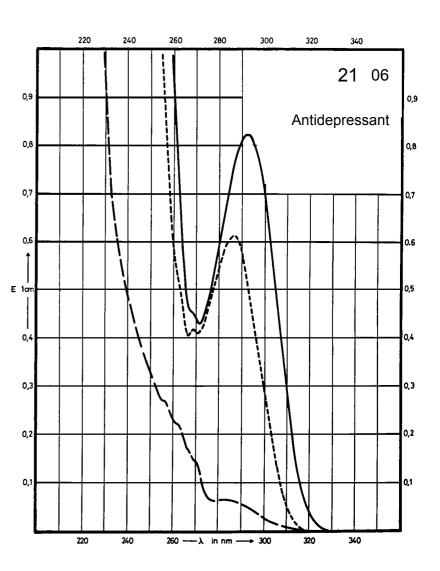


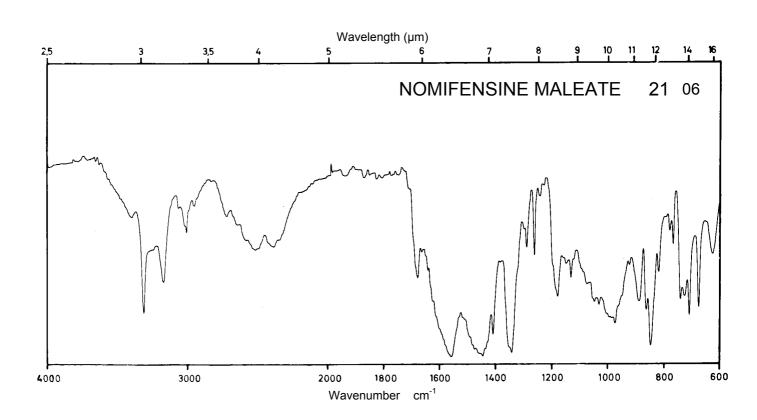

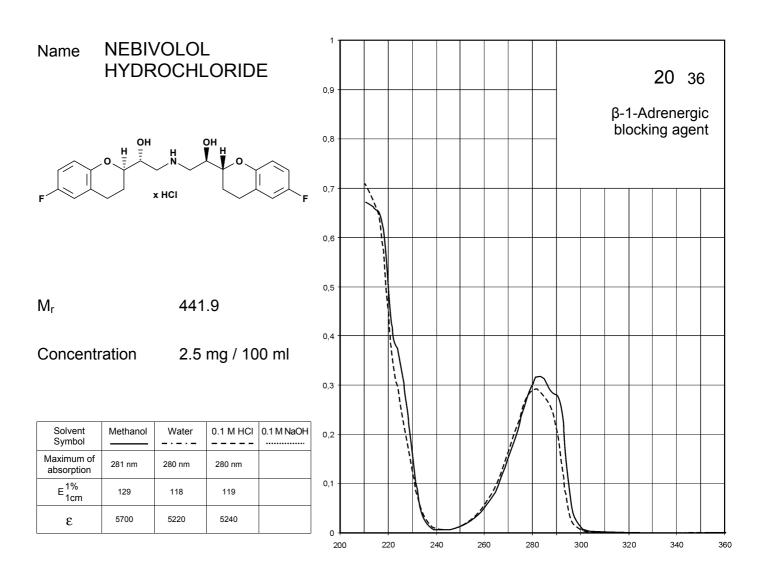

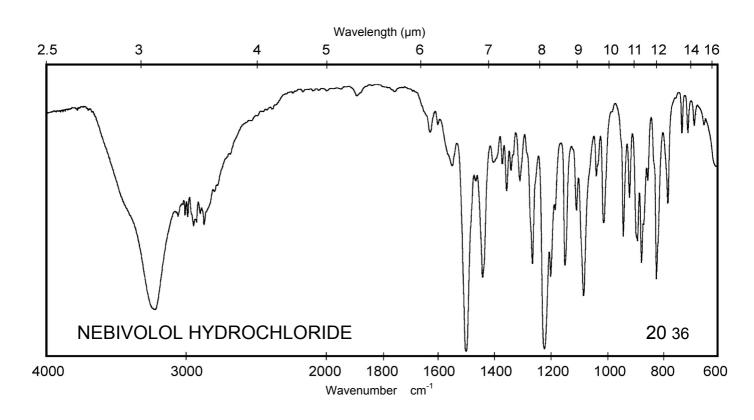





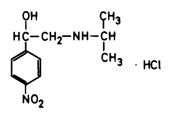


 M_r 260.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm		273 nm	278 nm
E ^{1%} 1cm	398		382	362
3	10370		9950	9430

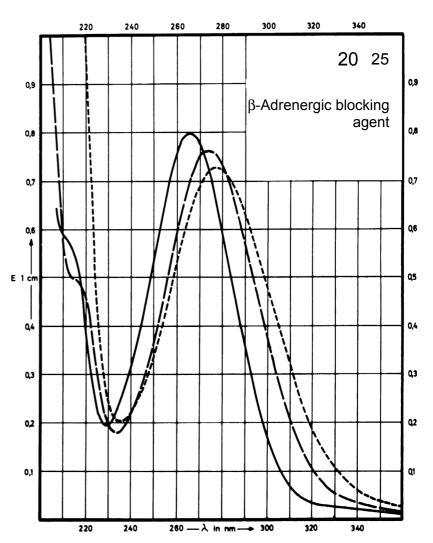

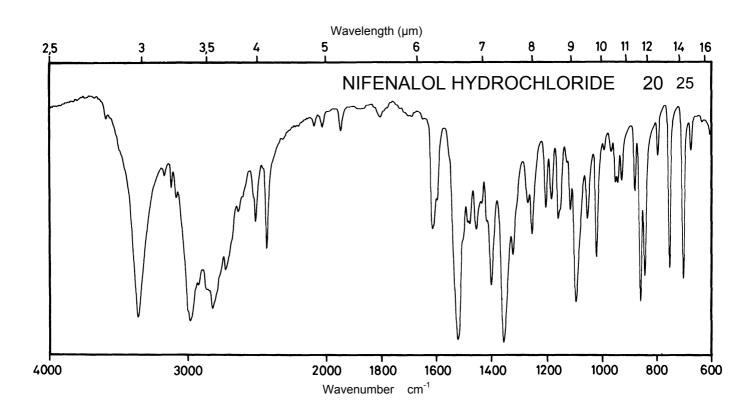




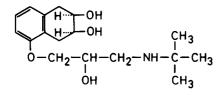


 M_{r} 309.4


Concentration 2 mg / 100 ml 10 mg / 100 ml

[
Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	279 nm 271 nm		277 nm 270 nm	278 nm 271 nm
E ^{1%} 1cm	36.6 35.5		32.3 32.8	33.2 34.1
ε	1130 1100		1000 1020	1030 1060

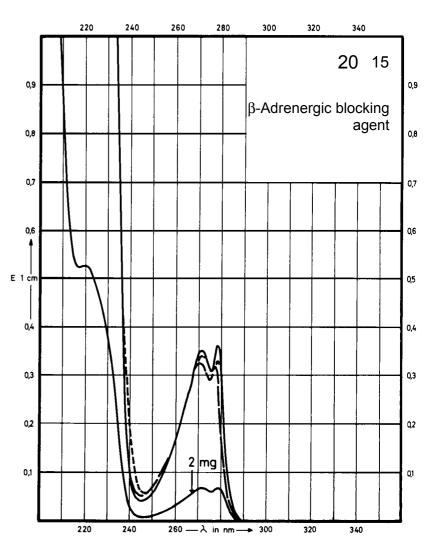

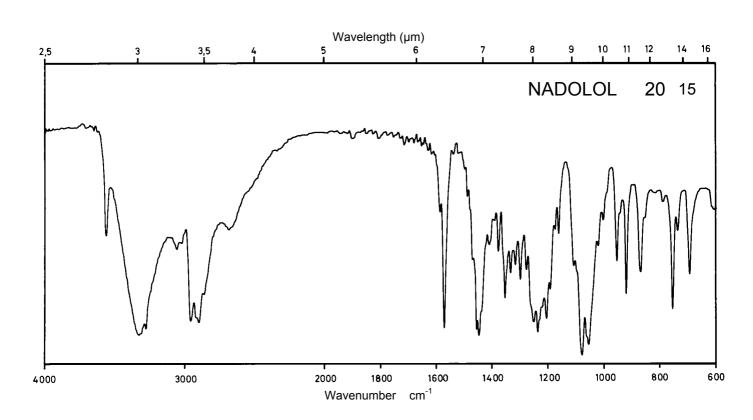





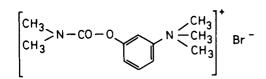

303.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm 261 nm		266 nm 260 nm	Decom- position observed
E ^{1%} 1cm	16.2 18.5		13.6 15.9	
3	490 560		410 480	

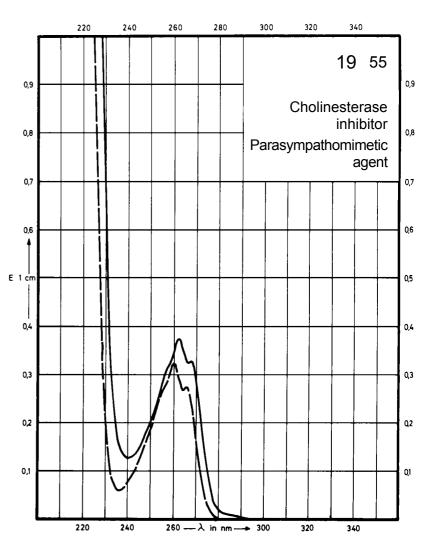


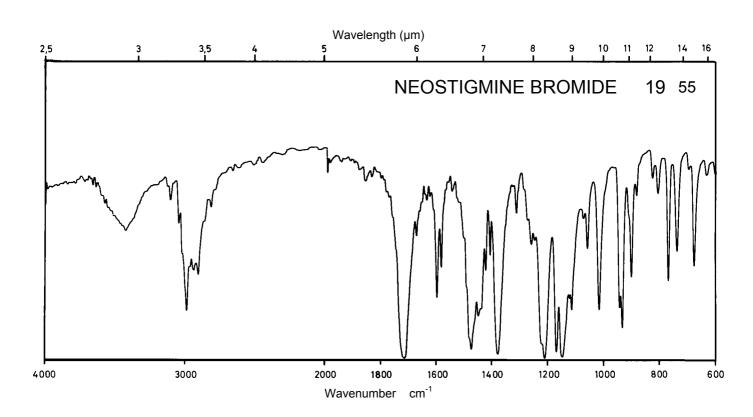




334.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm 266 nm		260 nm 266 nm	242 nm 295 nm
E ^{1%} 1cm	16 14		16 14	Decom- position observed
ε	535 470		535 470	

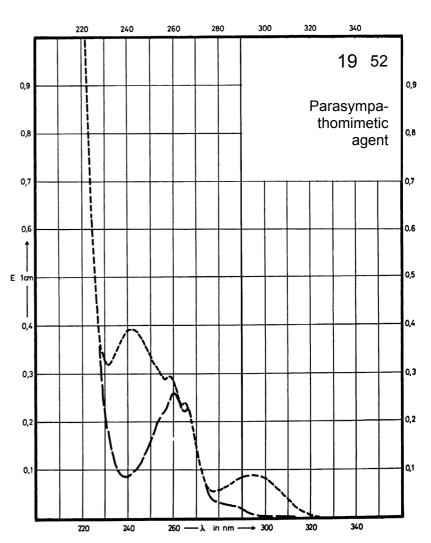

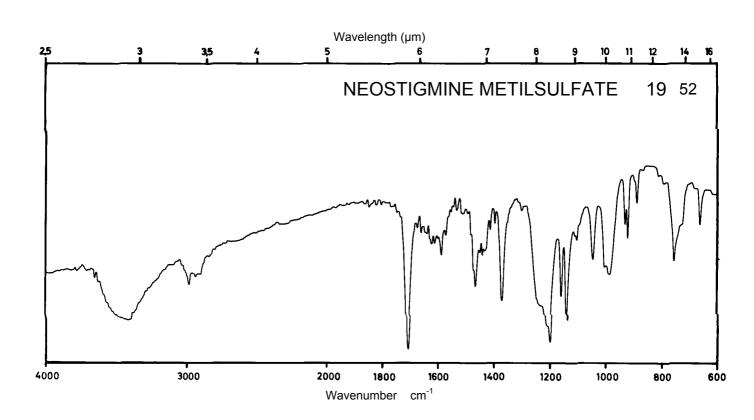








123.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm		260 nm	262 nm
E ^{1%} 1cm	225		420	260
ε	2770		5170	3200

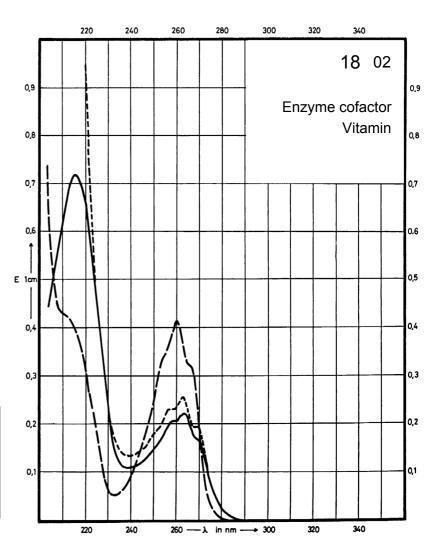

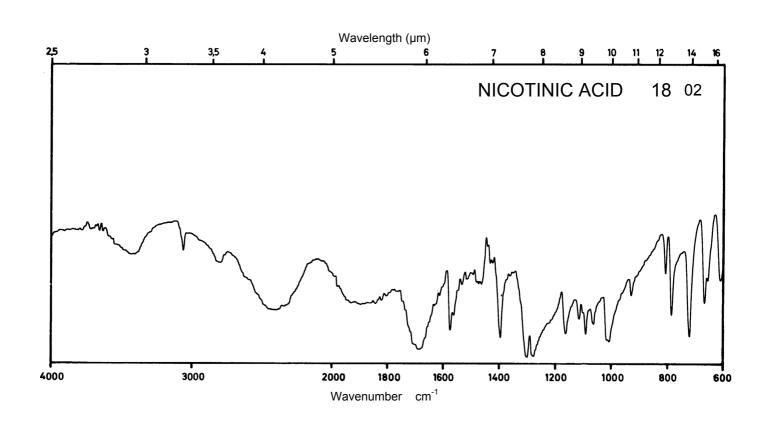





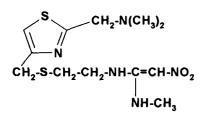


 M_r 331.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	326 nm 244 nm	314 nm 256 nm	316 nm	311 nm
E ^{1%} 1cm	255 578	475 241	36.7	389
3	19200 8400	15800 8000	1200	12900

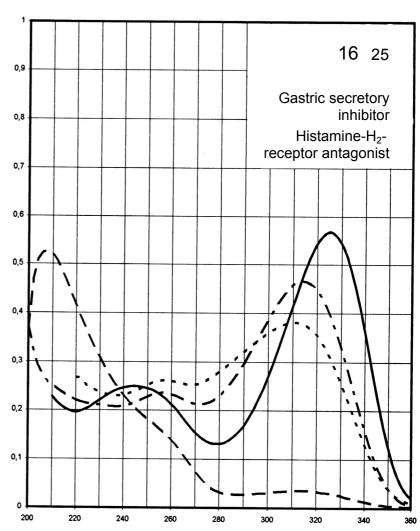

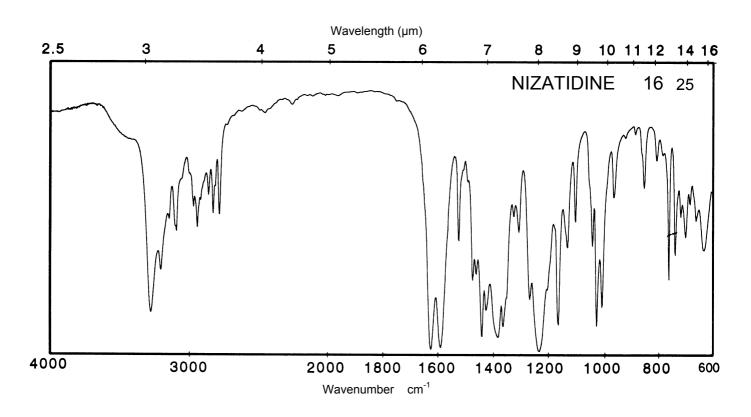




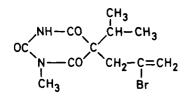


 M_r


303.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			222 nm	247 nm
E ^{1%} 1cm			235	243
3			7130	7350

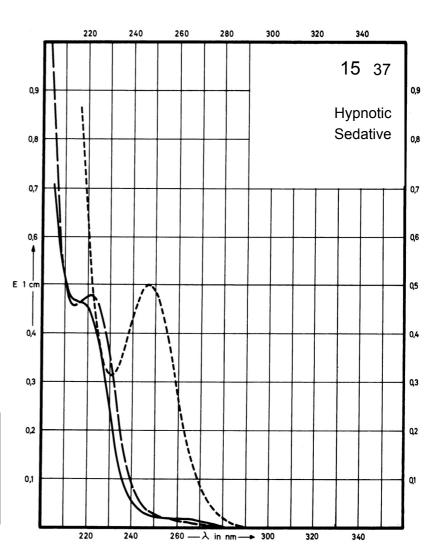

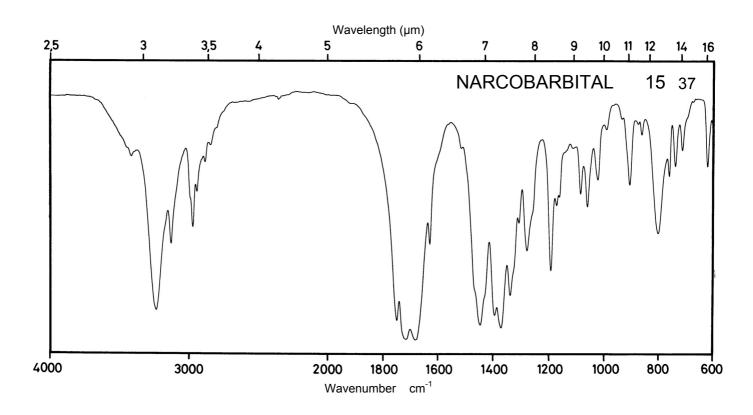

Name NITRAZEPAM

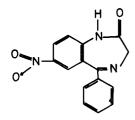


 M_{r}


281.3

Concentration 1


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	306 nm 257 nm		Decom- position observed	Decom- position observed
E ^{1%} 1cm	380 600			
ε	10690 16880			

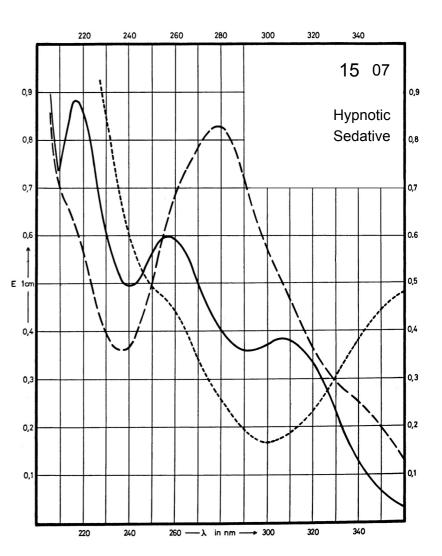

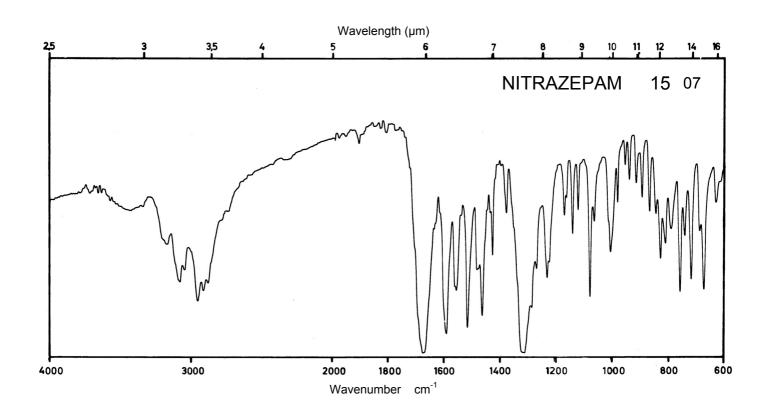





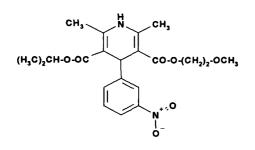


418.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	352 nm 231 nm	360 nm 236 nm		
E ^{1%} _{1cm}	174 815	161 545		
3	7270 34100	6700 22800		

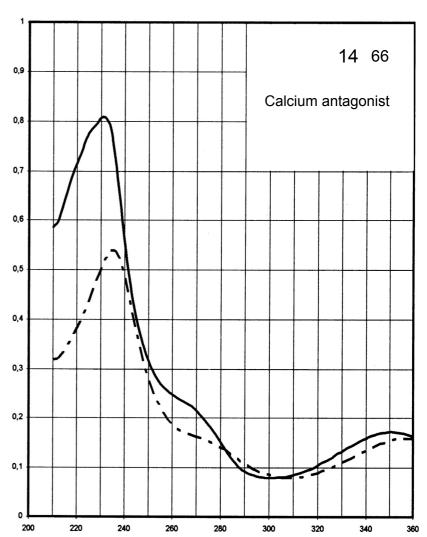

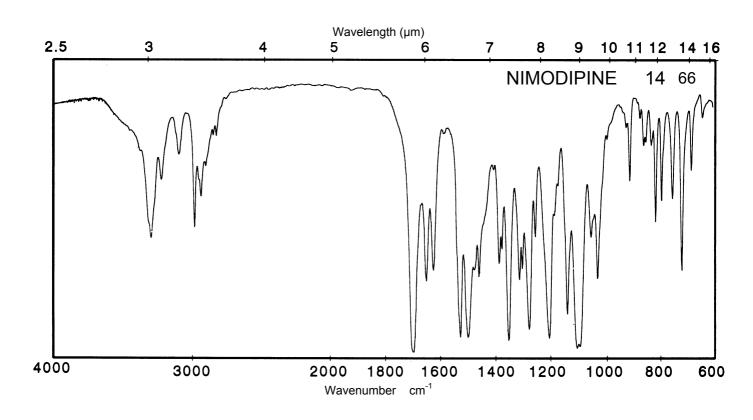




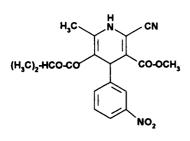


 M_r


385.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	374 nm 236 nm	381 nm 241 nm		
E ^{1%} 1cm	126 702	116 580		
3	4840 27050	4470 22350		

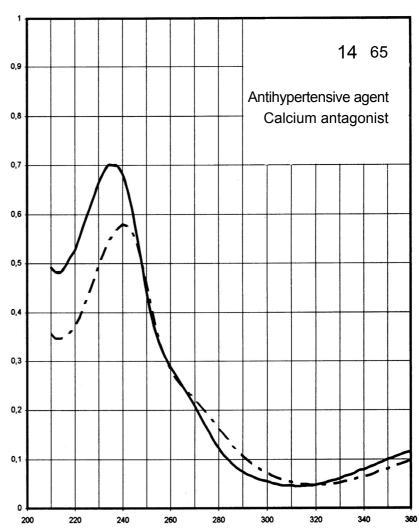

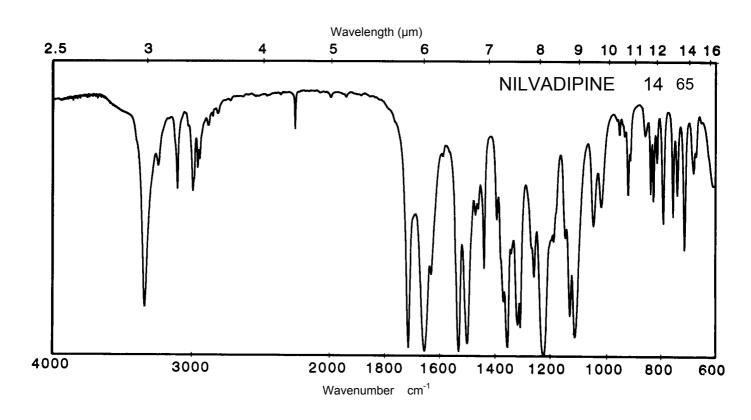




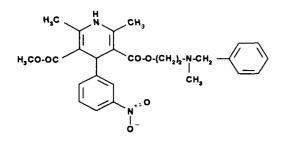


 M_{r}

516.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	355 nm	355 nm 239 nm	361 nm	
E ^{1%} 1cm	130	136 472	136	
3	6700	7000 24300	7050	

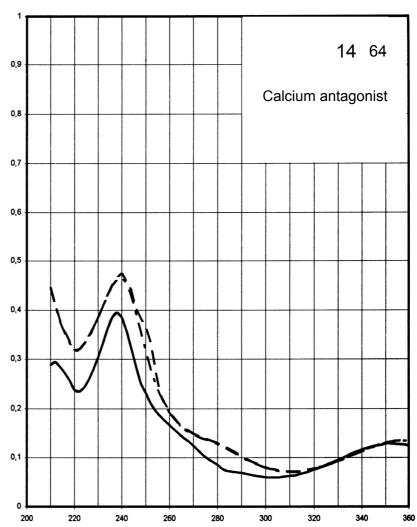

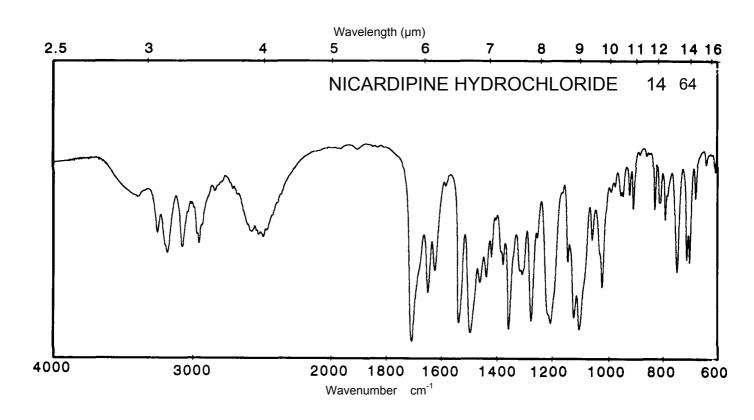



Name NITRENDIPINE


$$H_3CO-OC$$
 $CO-O-C_2H_5$
 NO_2


 M_r 360.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	352 nm 235 nm	360 nm 238 nm		
E ^{1%} 1cm	187 742	186 653		
3	6740 26700	6700 23500		

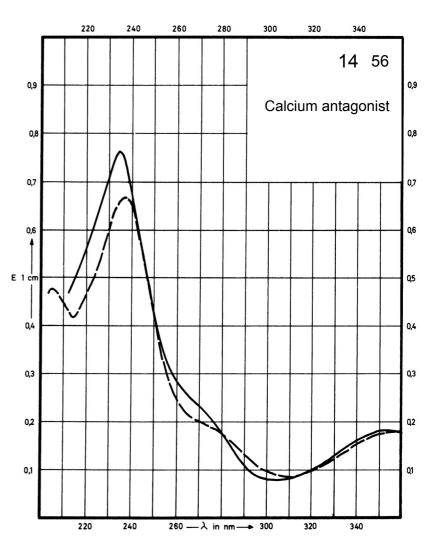

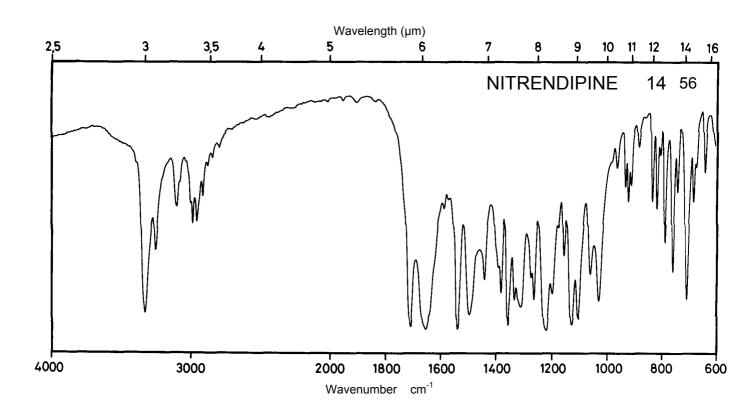

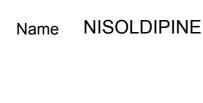

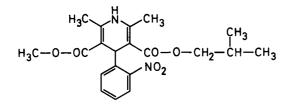

388.4

 M_r


Concentration 1 mg / 100 ml


 Solvent Symbol
 Methanol
 Water
 0.1 M HCI
 0.1 M NaOH

 Maximum of absorption
 345 nm 236 nm


 E 1% 1cm
 128 543

 E
 5000 21100

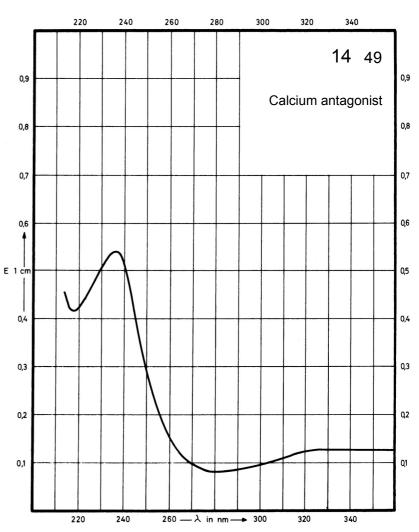

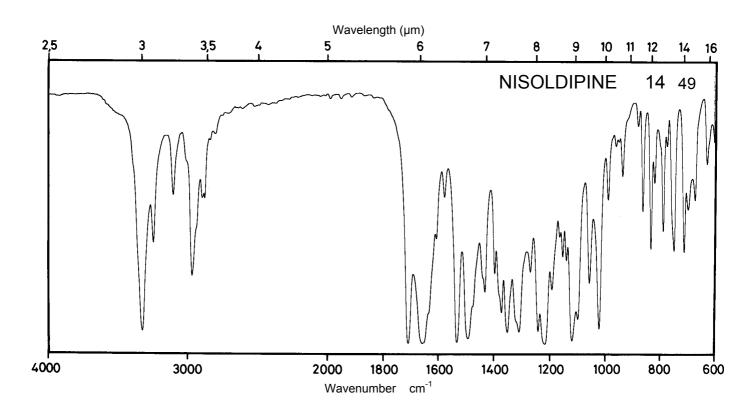




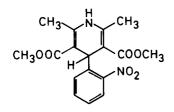


346.3


Concentration 1 mg / 100 ml 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	340 nm 235 nm		338 nm 238 nm	340 nm 238 nm
E ^{1%} 1cm	145 624		165 595	165 592
3	5010 21590		5740 20600	5740 20510

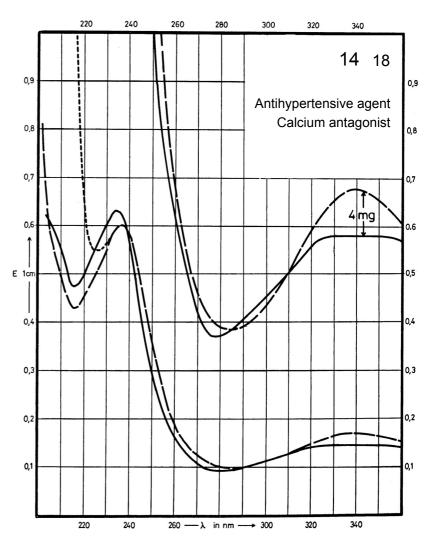

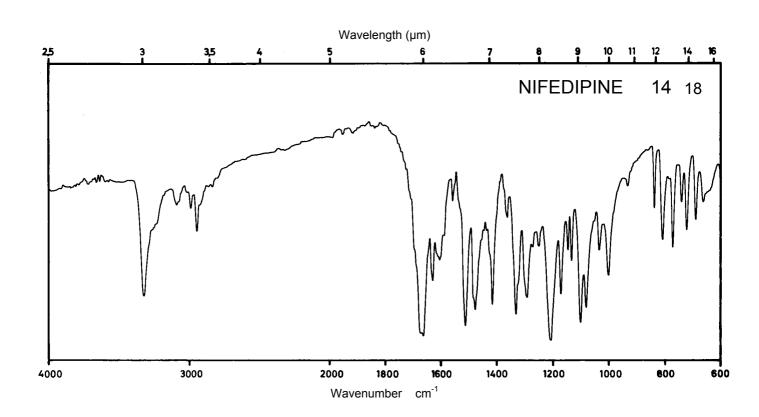



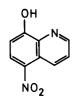



190.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	440 nm 355 nm 238 nm		370 nm 295 nm 248 nm	448 nm 328 nm 247 nm
E ^{1%} 1cm	214 439 1238		386 283 1685	1621 153 809
ε	4060 8330 23540		7350 5380 30050	30830 2910 15390

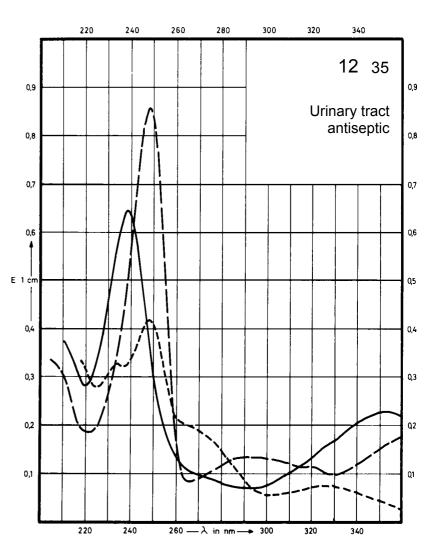

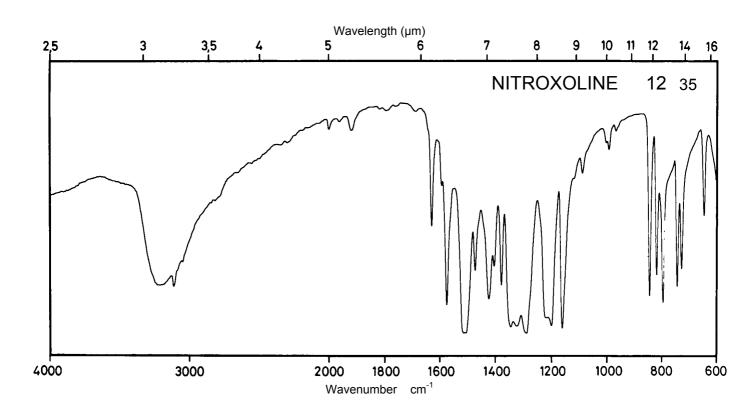



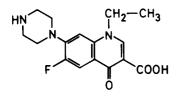



319.3

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			315 nm 277 nm	335 nm 323 nm 272 nm
E 1% 1cm			383 1350	430 440 1130
ε			12200 43100	13700 14000 36100

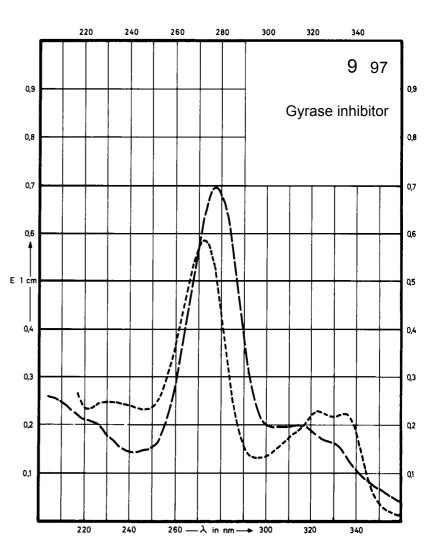

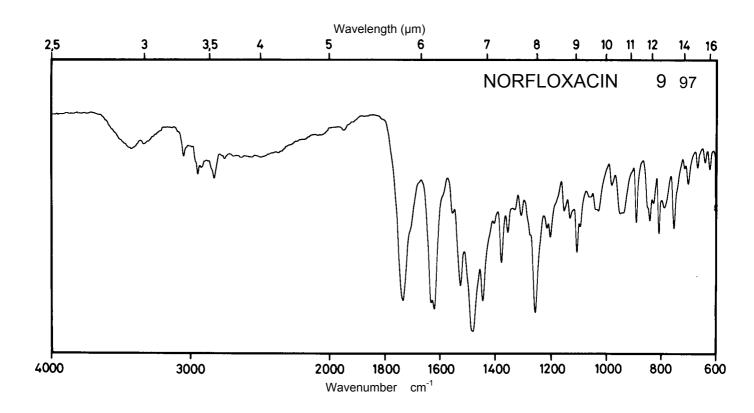



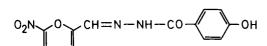



 M_r


275.2

Concentration


1 mg / 100 ml

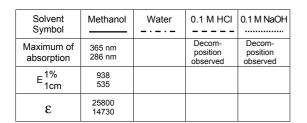

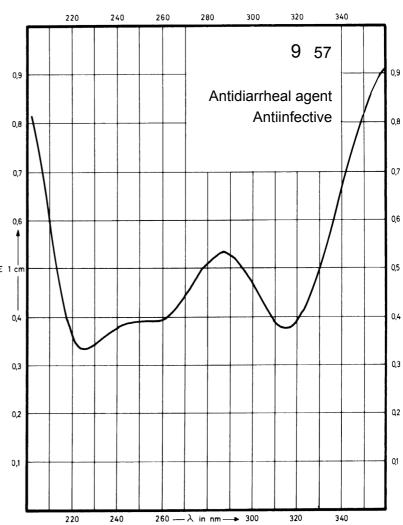

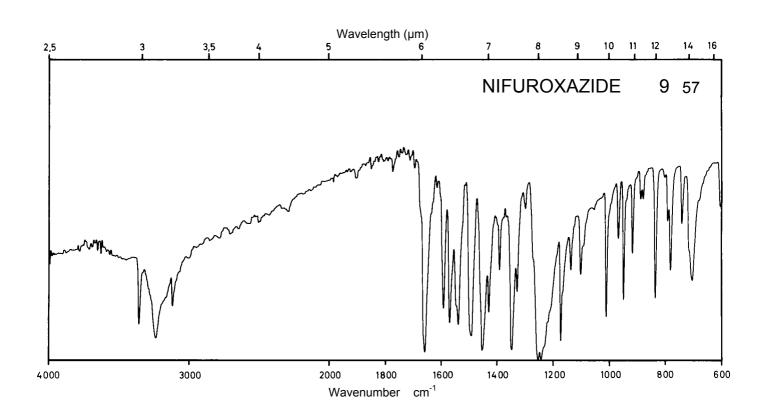





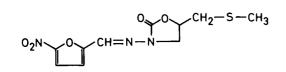

 M_r 285.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	355 nm 262 nm		367 nm 259 nm	
E ^{1%} 1cm	657 348		615 453	
3	18760 9910		17540 12930	

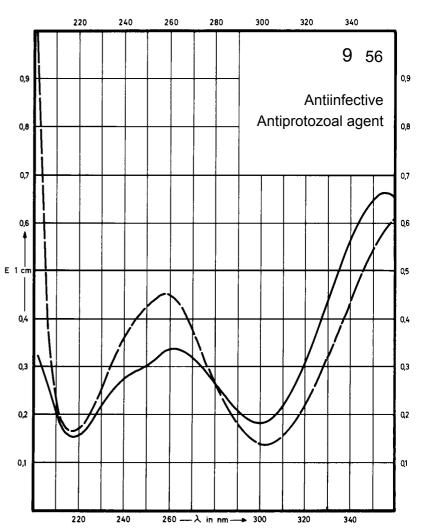

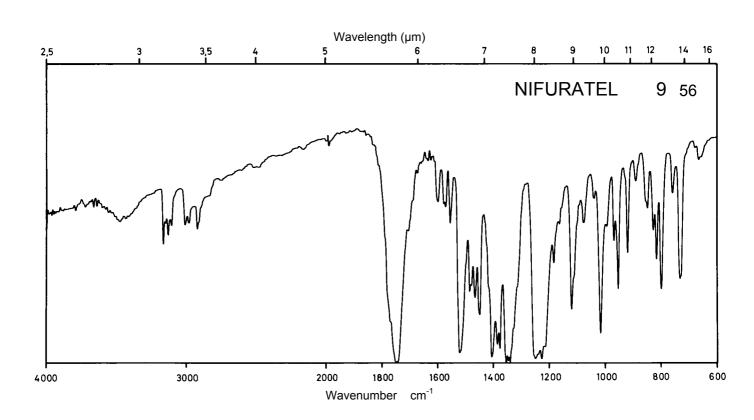

Name NALIDIXIC ACID

СН₃ СН₂ СООН


 M_{r}


232.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	329 nm 258 nm		315 nm 257 nm	334 nm 258 nm
E ^{1%} 1cm	496 1110		480 1233	495 1142
ε	11520 25780		11230 28630	11490 26510

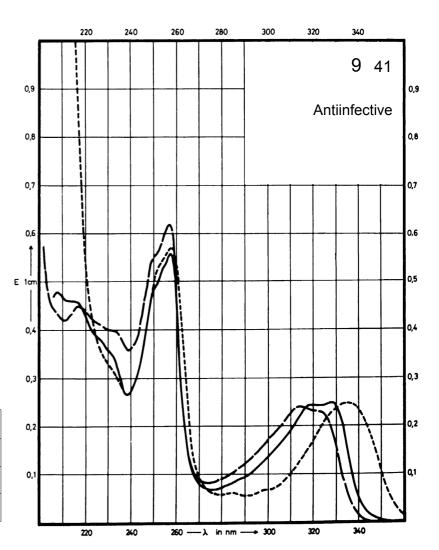

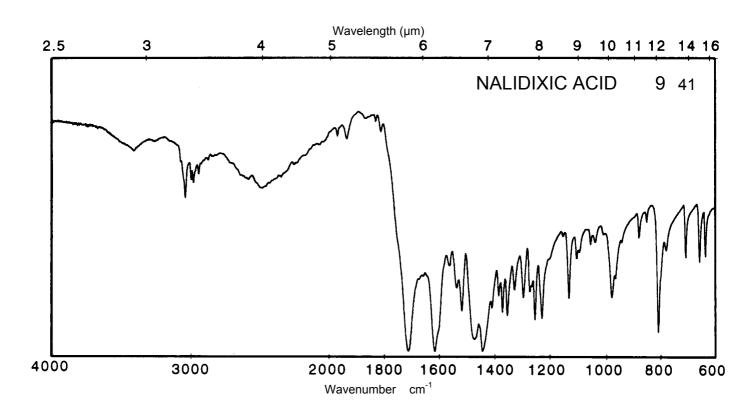


Name NITROFURANTOIN

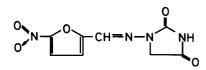


 M_{r}


238.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	357 nm 269 nm		266 nm	281 nm
E ^{1%} _{1cm}	669 427		567	568
ε	15940 10170		13510	13530

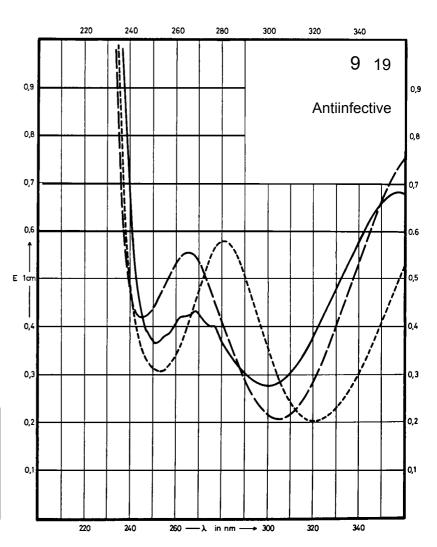

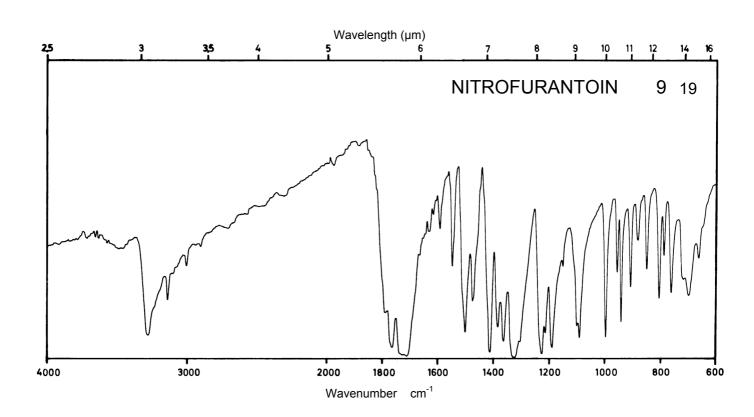


Name NITROFURAL

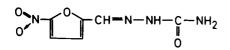


 M_{r}


198.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm		261 nm	284 nm
E ^{1%} 1cm	563		610	475
3	11160		12090	9410

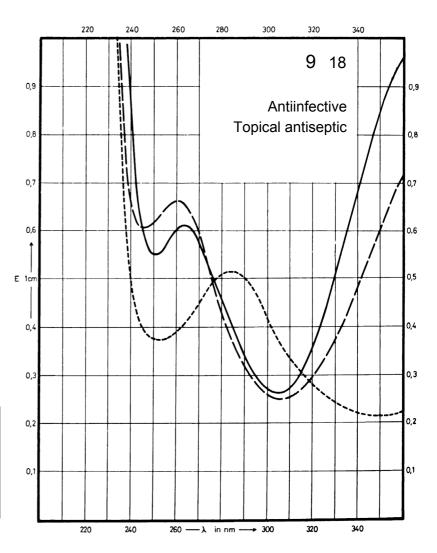

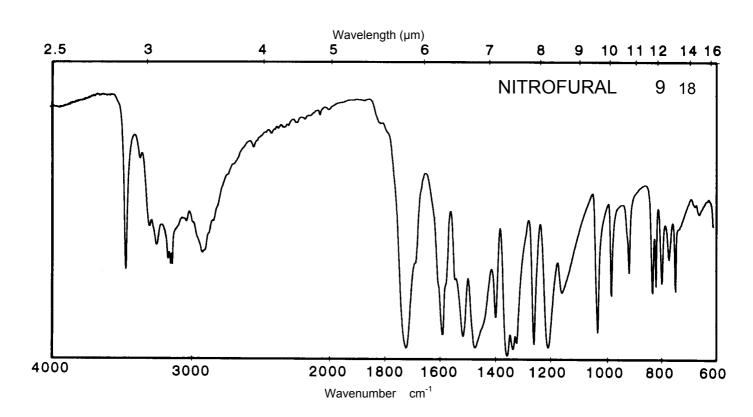



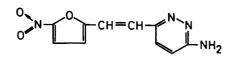



232.2

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	374 nm 288 nm 253 nm		383 nm 284 nm 262 nm	398 nm 294 nm 248 nm
E ^{1%} 1cm	908 610 446		914 624 582	689 601 394
ε	21080 14160 10360		21220 14490 13510	16000 13960 9150

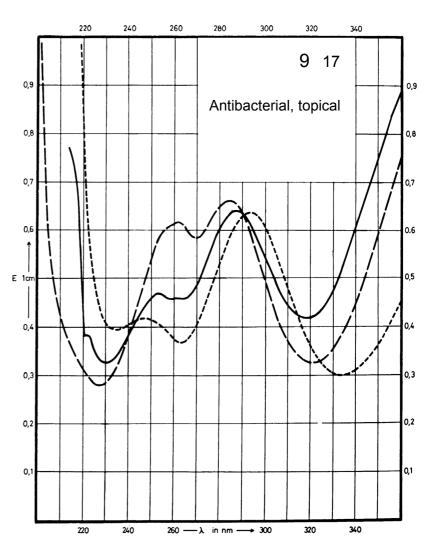

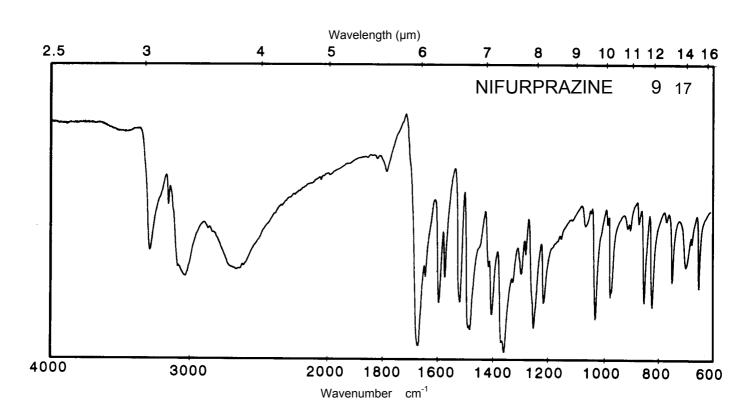



Name NOSCAPINE HYDROCHLORIDE

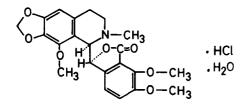


 M_{r} 467.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm		312 nm	280 nm
E ^{1%} 1cm	94		83	72
ε	4400		3880	3370

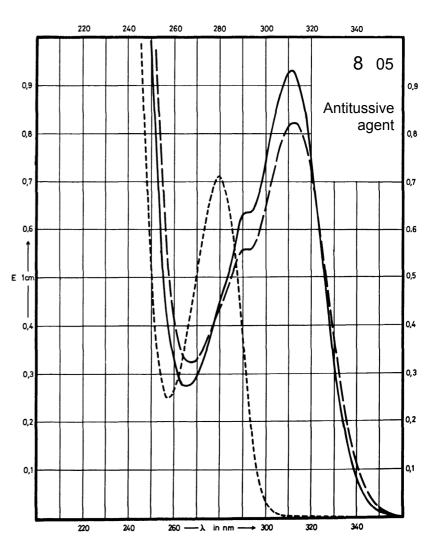

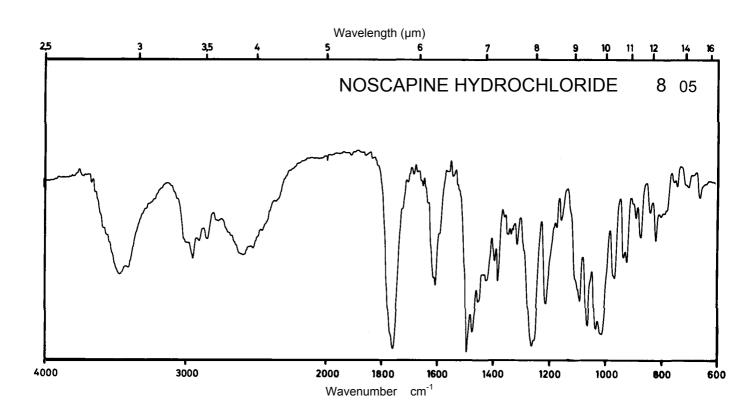





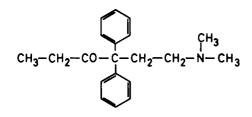


 M_r 295.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm 260 nm			
E ^{1%} 1cm	13.8 14.7			
ε	410 435			

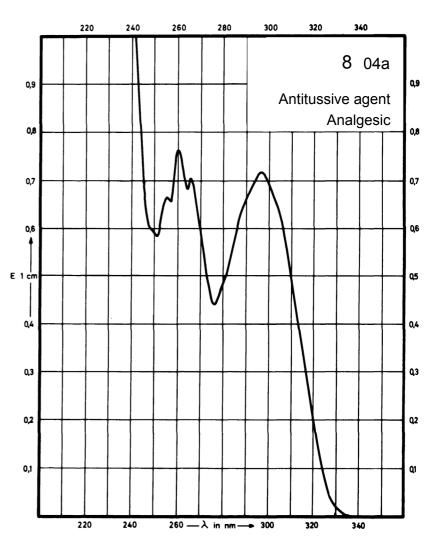

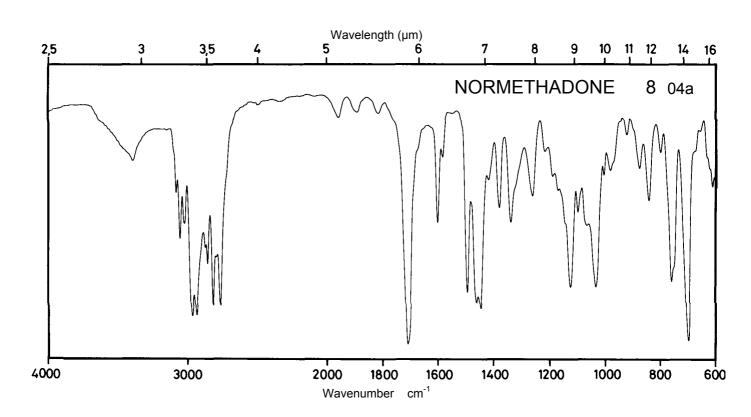



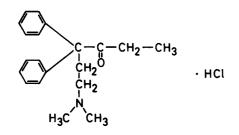

 M_r


331.9


Concentration 27 mg / 100 ml (CH₃ŎH)

33 mg /100 ml


 (H_2O)

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	294 nm 258 nm	292 nm 258 nm		
E ^{1%} 1cm	13.3 13.5	16.6 16.0		
ε	440 450	550 530		

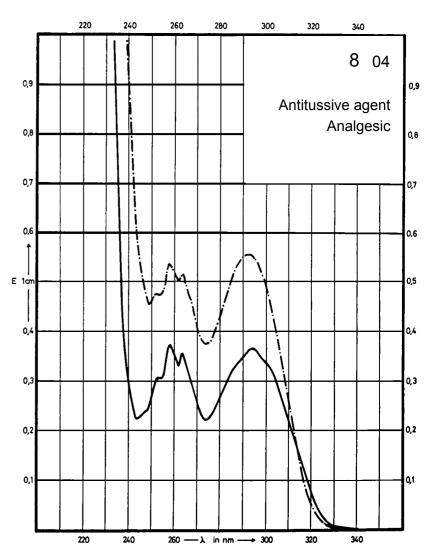

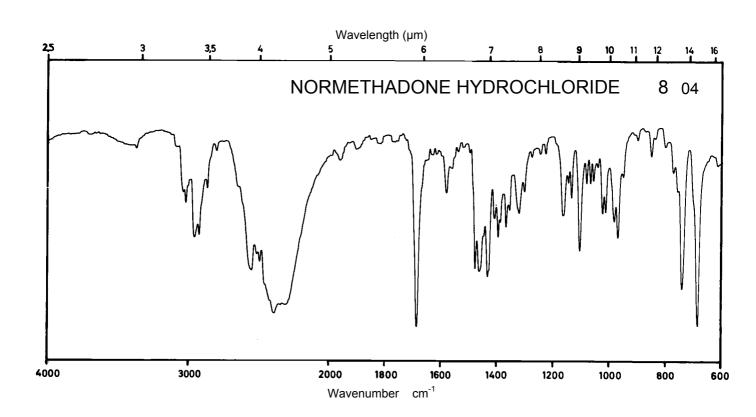




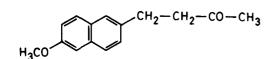


228.3


Concentration 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	333 nm 270 nm 260 nm			
E ^{1%} 1cm	88 216 219			
ε	2000 4940 5000			

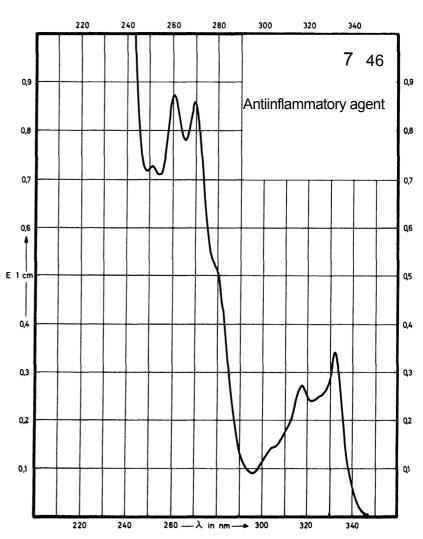

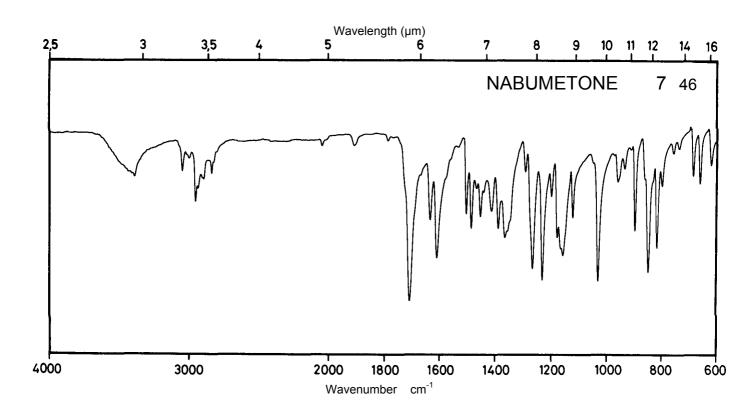




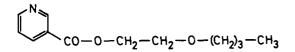


223.3


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm		260 nm	
E ^{1%} 1cm	134		249	
ε	2970		5560	

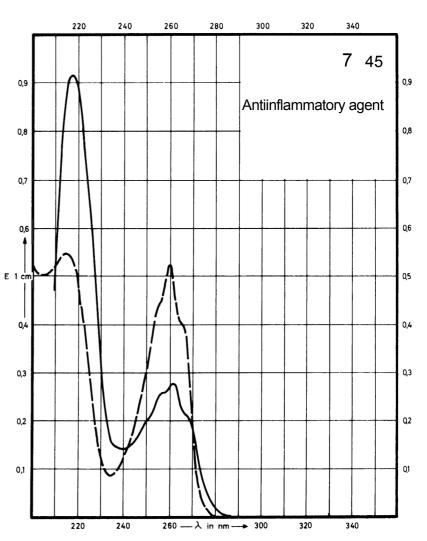

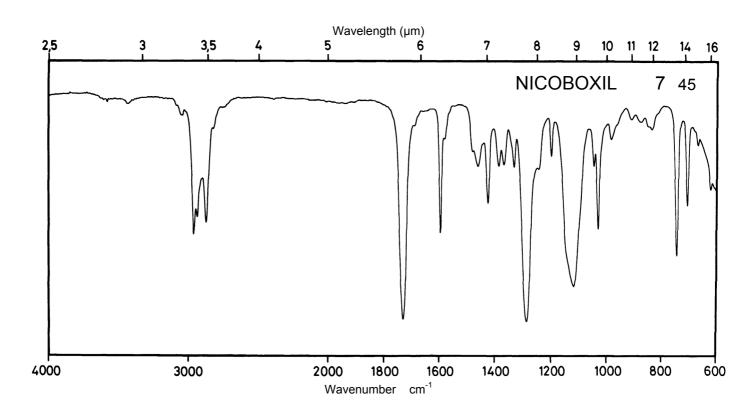





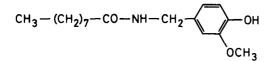


 M_{r} 293.4

Concentration 2 mg / 100 ml 3 mg / 100 ml

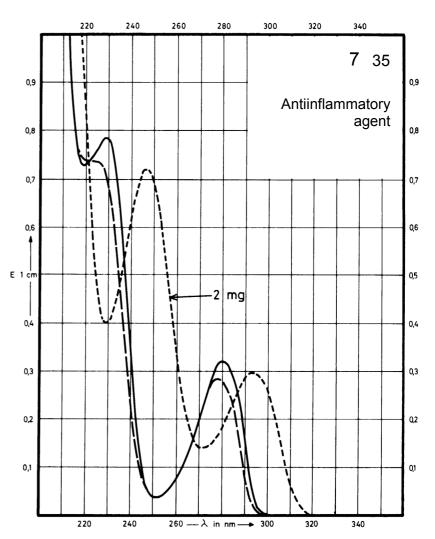

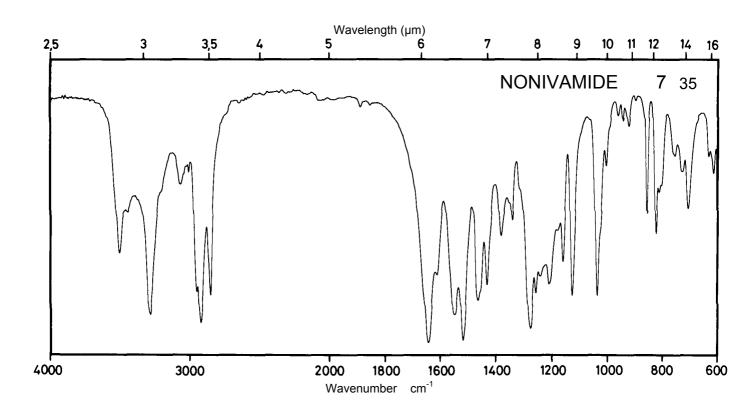




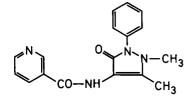


 M_r 308.3


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm		260 nm	
E ^{1%} 1cm	442		503	
3	13600		15500	

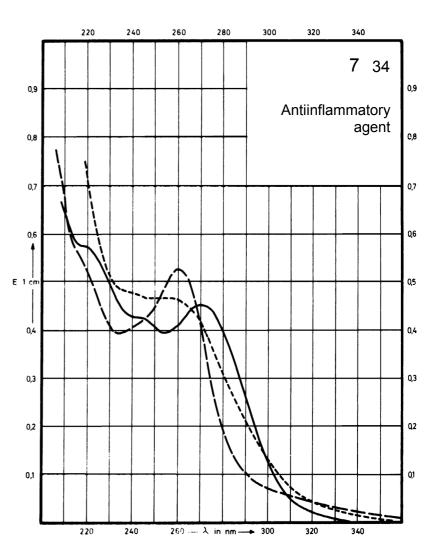

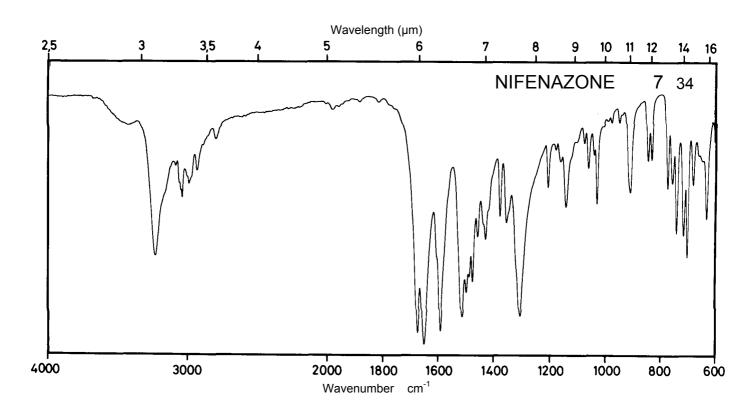



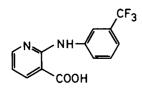



Concentration

282.2


0.5 mg / 100 ml 3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	342 nm 289 nm		329 nm 255 nm	288 nm
E ^{1%} 1cm	192 987		202 471	788
3	5410 27870		5690 13300	22250

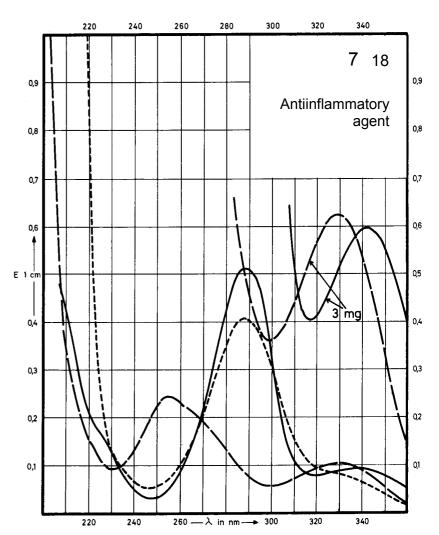

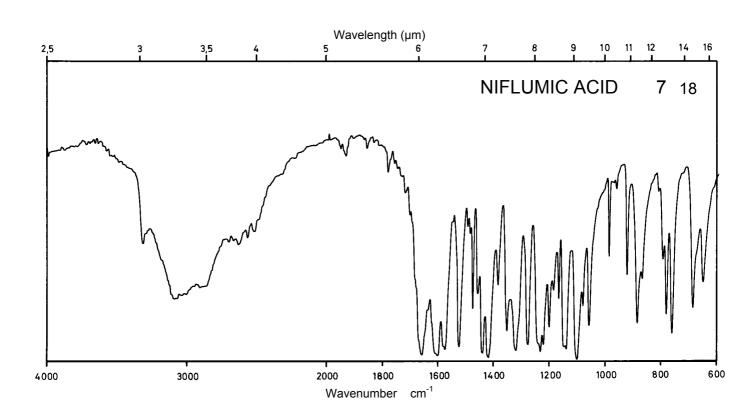




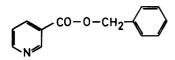


213.4


Concentration 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm			
E ^{1%} 1cm	152			
ε	3140			

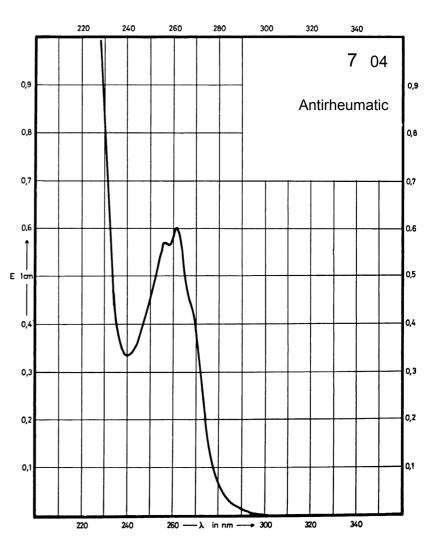

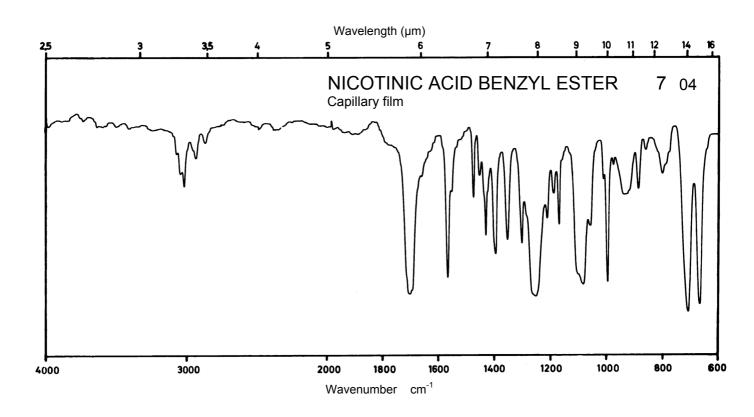



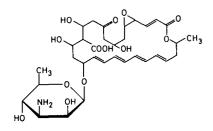

 M_r 665.8


Concentration 0.5 mg / 100 ml (846 μg / mg)


 Solvent Symbol
 Methanol ————
 Water ————
 0.1 M HCI ————
 0.1 M NaOH


 Maximum of absorption
 317 nm 302 nm 289 nm
 319 nm 304 nm 291 nm


 E 1% 1cm
 1187 1270 815
 1070 1173 758

 T 79000 84500
 71300 78100 54300
 78100 50500

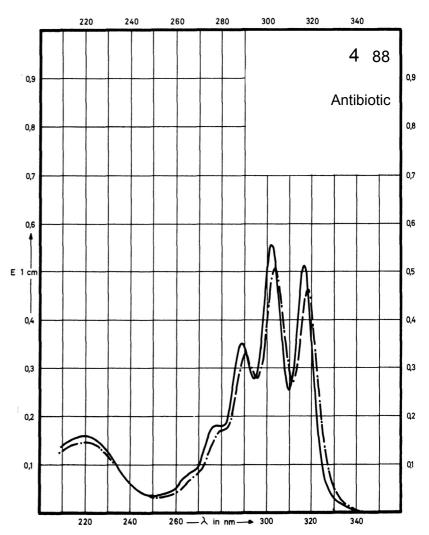


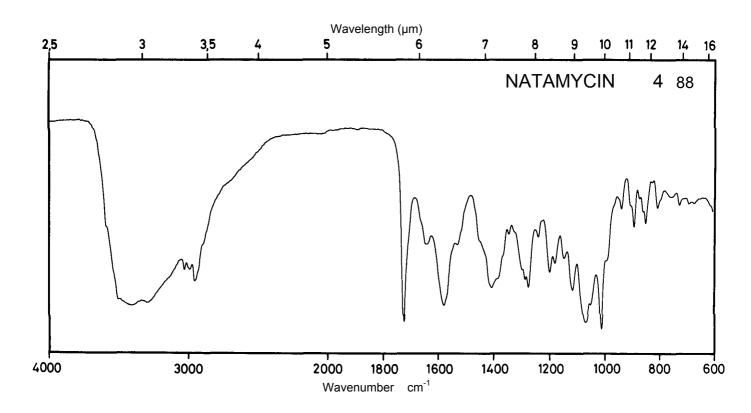

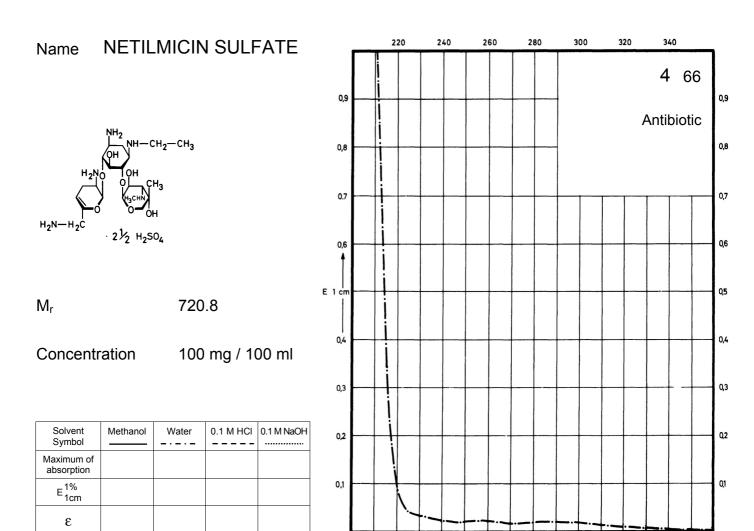

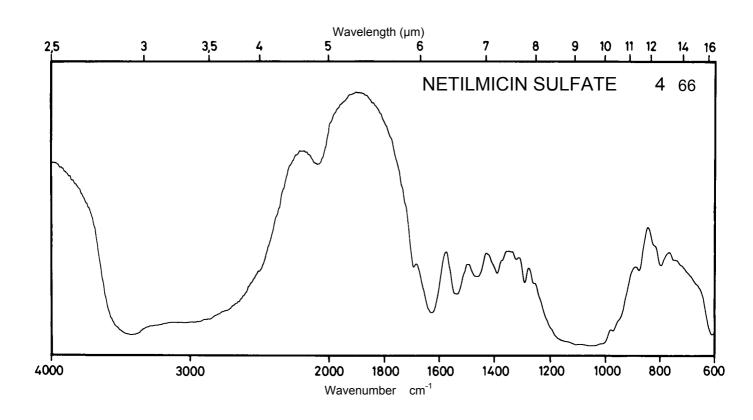



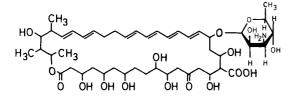

 M_r


926.1


1 mg / 100 ml Concentration

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 320 nm 305 nm 231 nm 475 530 320 nm 305 nm 292 nm 318 nm 303 nm 230 nm Maximum of absorption 600 660 530 600 E^{1%} 1cm 230 43990 405 49080 55570 37510 ε 61120 49080 21300


Further maxima (E1%,1 cm):
Methanol: 405 nm (22), 381 nm (21), 361 nm (15)
0.1 M HCl: 409 nm (15), 385 nm (16), 363 nm (13)
0.1 M NaOH: 407 nm (18), 384 nm (20), 363 nm (17)

22230

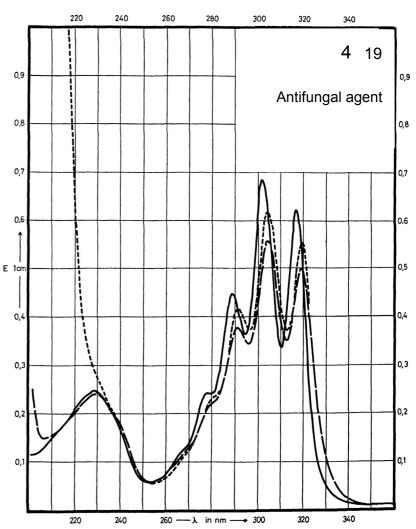

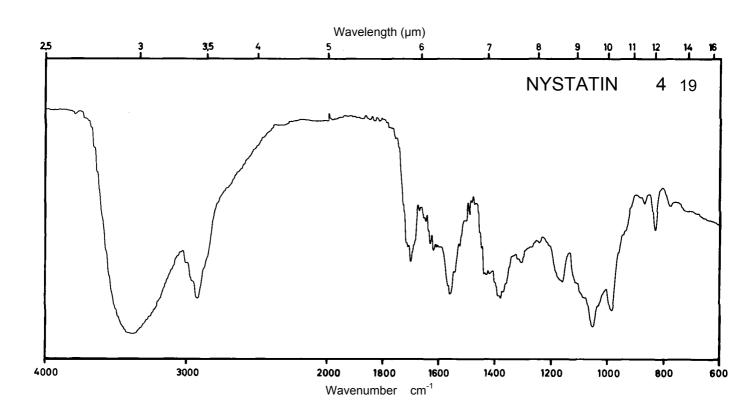



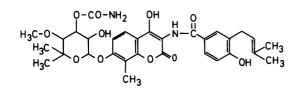



 M_r

612.6


Concentration 1 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	305 nm 248 nm			305 nm
E 1% 1cm	304 339			552
3	18620 20770			33820

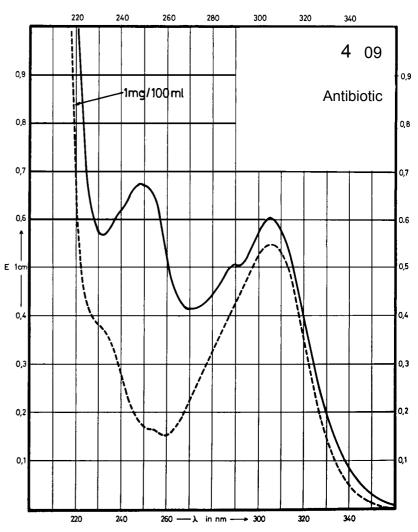

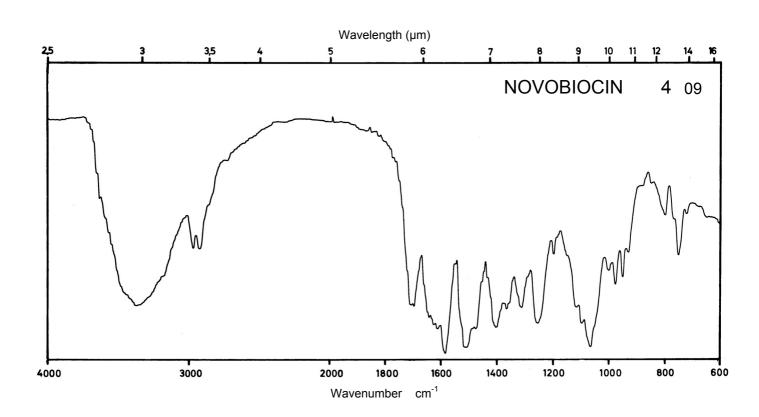



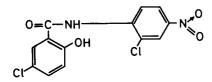



327.1

Concentration 1 r


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	332 nm			377 nm 334 nm
E ^{1%} _{1cm}	527			422 412
3	17240			13800 13480

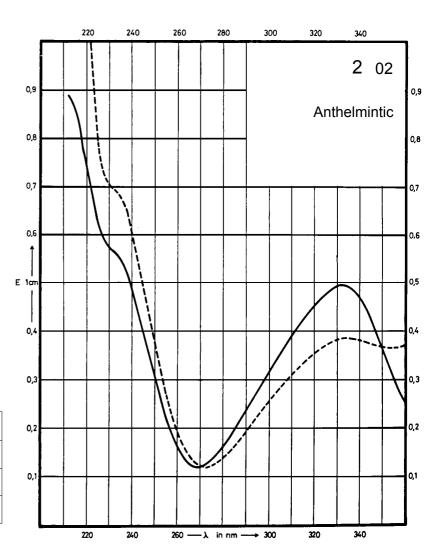

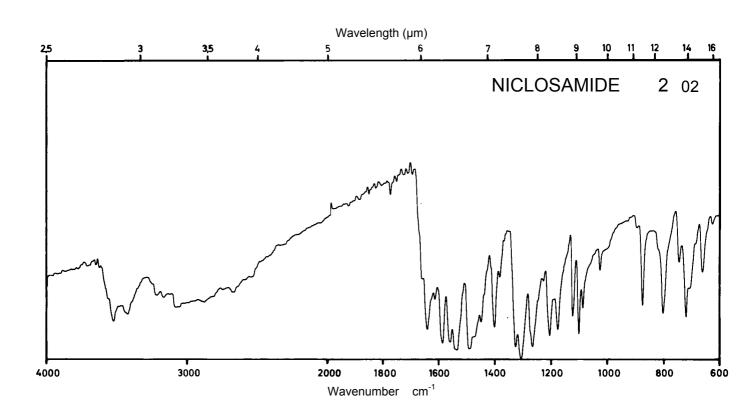

Name NEFOPAM HYDROCHLORIDE

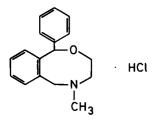



 M_r 289.8

Concentration 32 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI*	0.1 M NaOH*
Maximum of absorption	275 nm 267 nm		275 nm 267 nm	259 nm
E ^{1%} _{1cm}	20.4 24.4		21.6 25.5	18.0
ε	590 706		626 740	522

^{* 1} M HCl + Methanol (1+9) 1 M NaOH + Methanol (1+9)

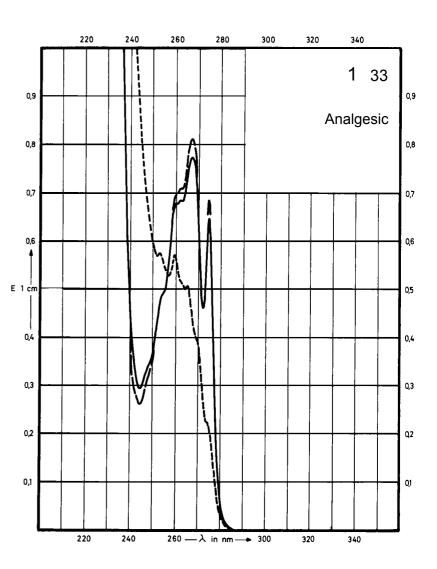

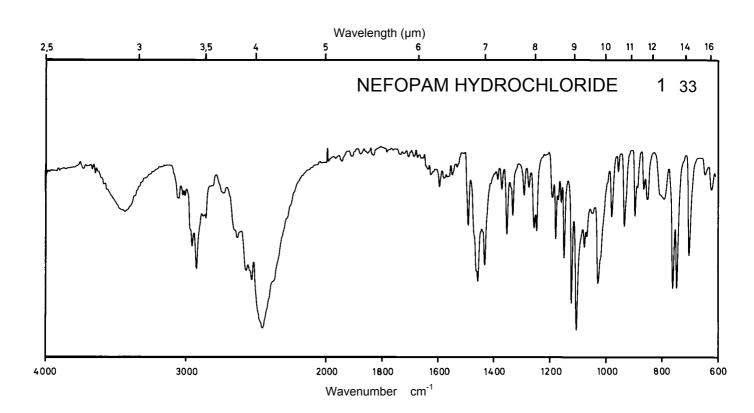



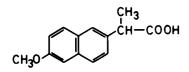



230.3

Concentration


2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	331 nm 316 nm 271; 262 nm		328 nm 315 nm 272; 262 nm	330 nm 316 nm 271; 261 nm
E 1% 1cm	80 63 225; 220		63 52 215; 208	70 52 218; 218
ε	1840 1450 5180; 5070		1450 1200 4950; 4790	1610 1200 5020; 5020

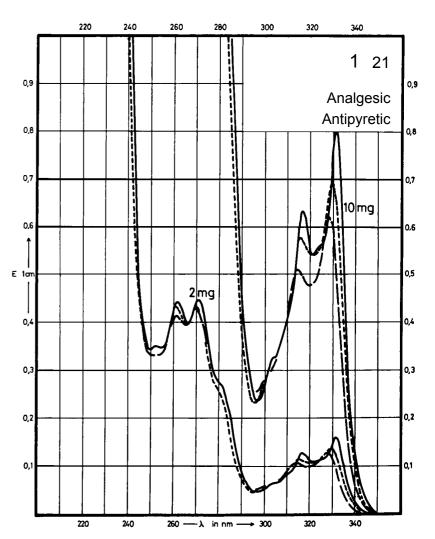

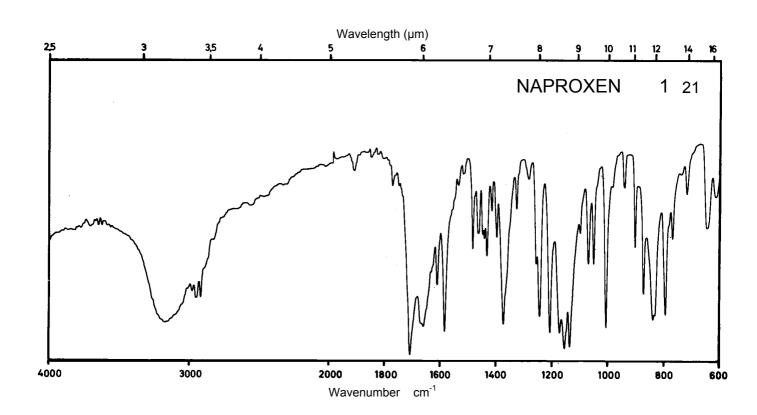



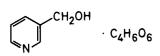



 M_r

259.2


Concentration 3 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm		260 nm	260 nm
E ^{1%} 1cm	113		206	112
3	2930		5350	2900

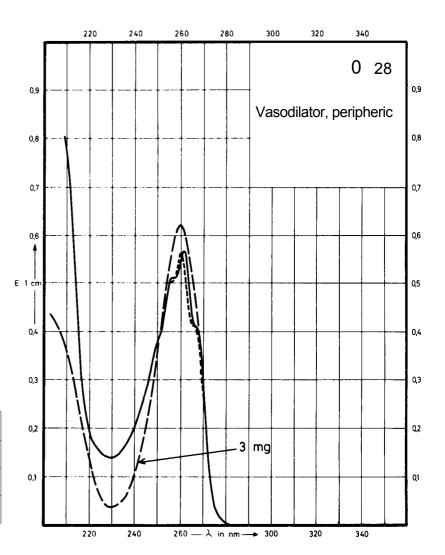

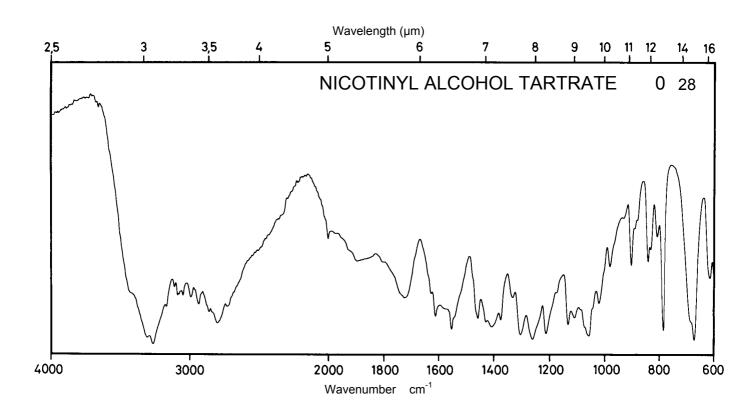




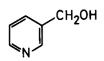


109.1


Concentration 1.3 mg / 100 ml 2.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm		260 nm	260 nm
E ^{1%} _{1cm}	264		509	269
3	2880		5550	2930

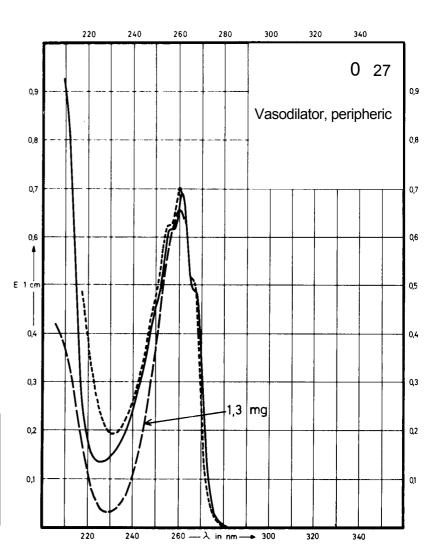

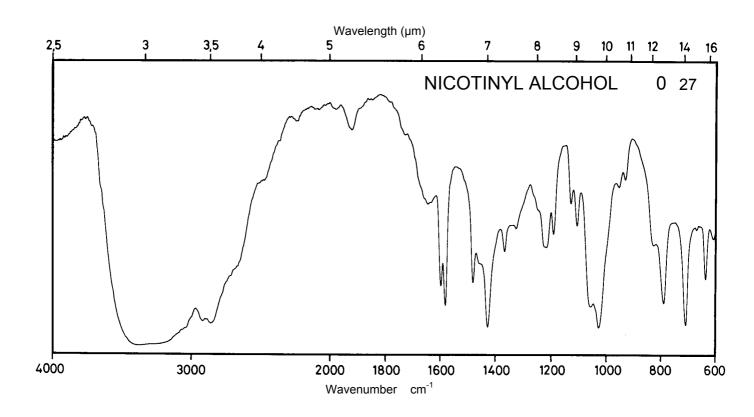




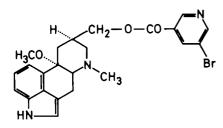


 M_r 484.4


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm		288 nm	283 nm
E ^{1%} 1cm	183		197	164
3	8800		9500	7900

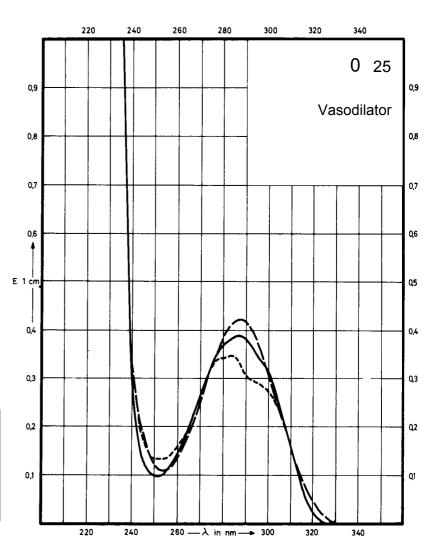

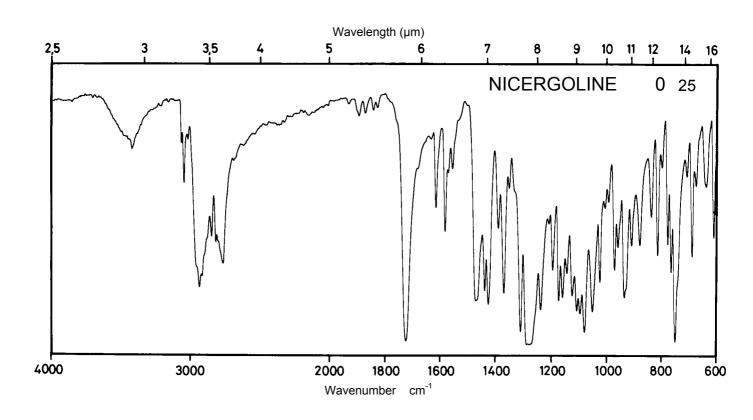




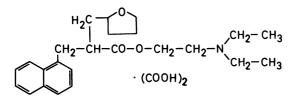


 M_r 473.6


Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm	283 nm	283 nm	283 nm
	273 nm	273 nm	273 nm	273 nm
E ^{1%}	155	144	144	142
1cm	128	121	121	119
ε	7340	6820	6820	6730
	6060	5730	5730	5640

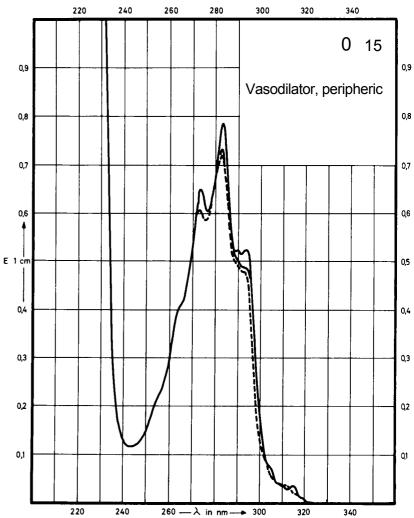

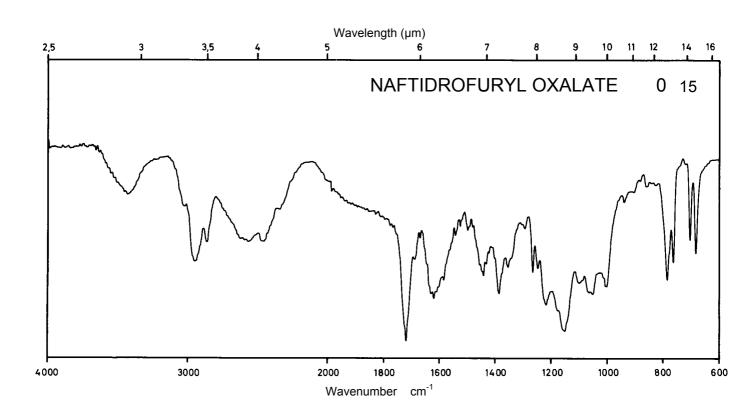




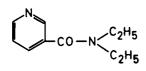


178.2


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm		264 nm	262 nm
E ^{1%} 1cm	184		293	192
ε	3280		5220	3420

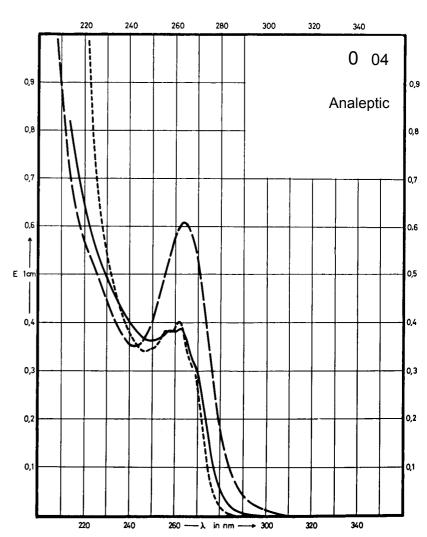

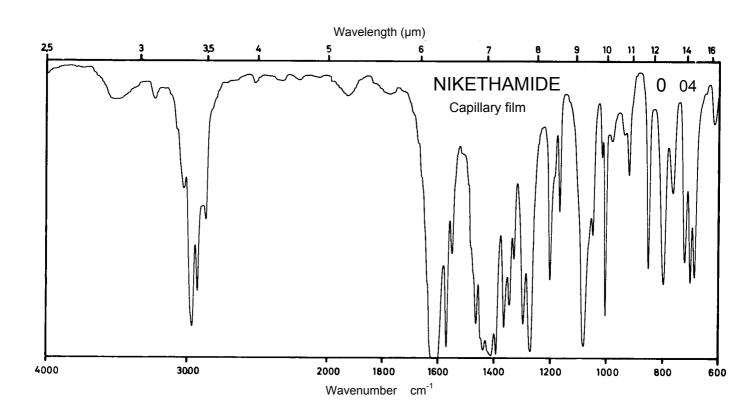

Name ORAZAMIDE

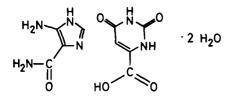



 M_{r}

318.3


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			271 nm 239 nm	279 nm
E ^{1%} 1cm			495 296	580
ε			15760 9400	18500

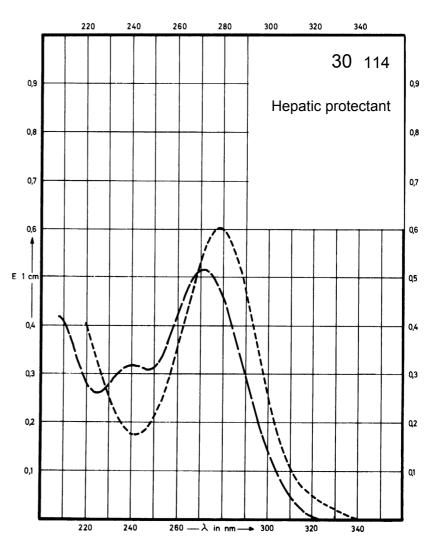

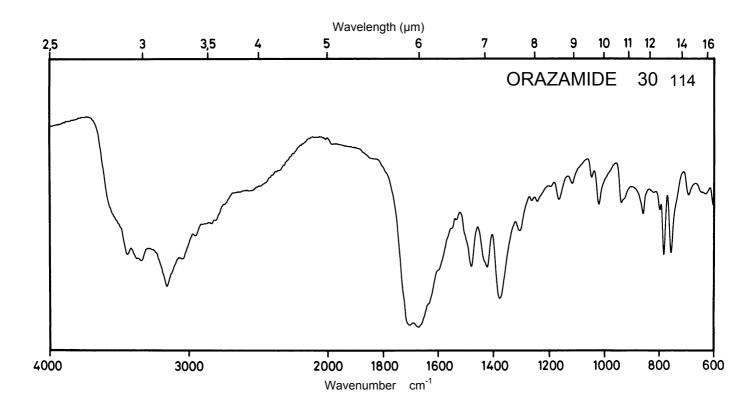




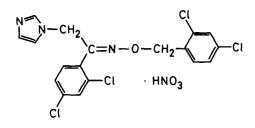


492.2


Concentration 1.2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

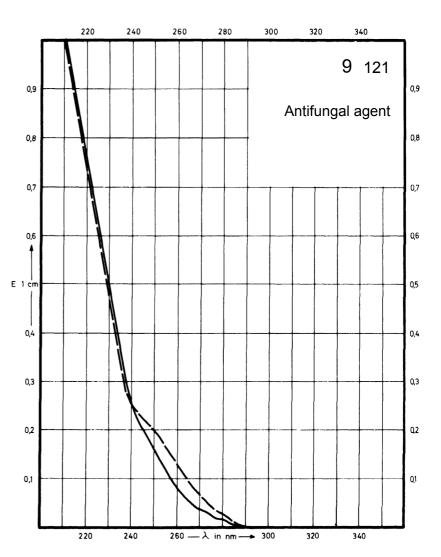

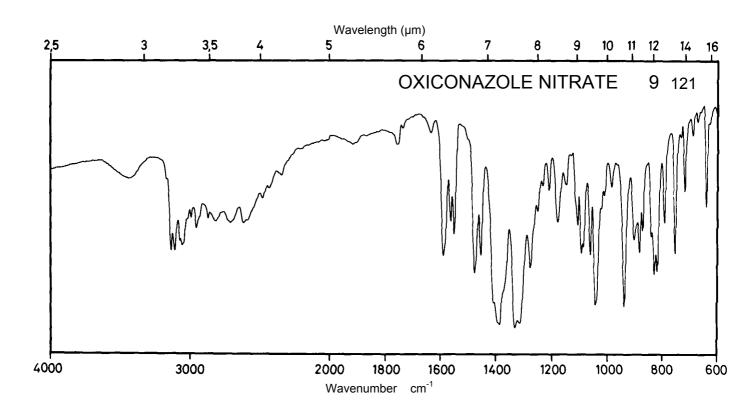




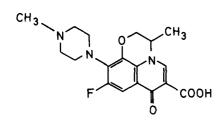


 M_r 361.4

Concentration 0.5 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	298 nm		294 nm	288 nm
E ^{1%} 1cm	995		910	740
3	36000		33000	26700

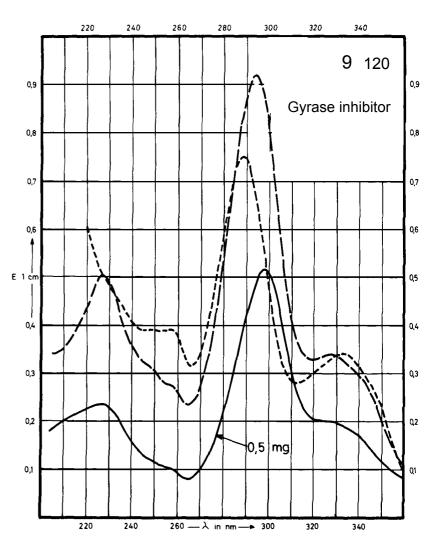

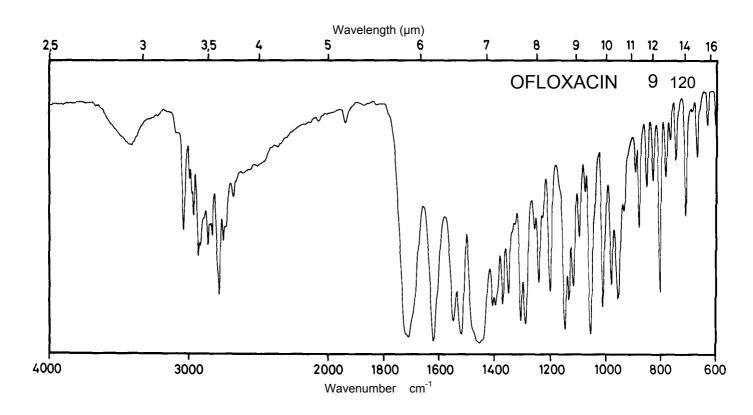

Name OXETORONE FUMARATE


 M_{r}


435.5

Concentration 1.

1.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	297 nm 254 nm	297 nm 252 nm	297 nm 254 nm	
E ^{1%} 1cm	458 278	443 276	443 260	
ε	19950 12100	19290 12020	19290 11290	

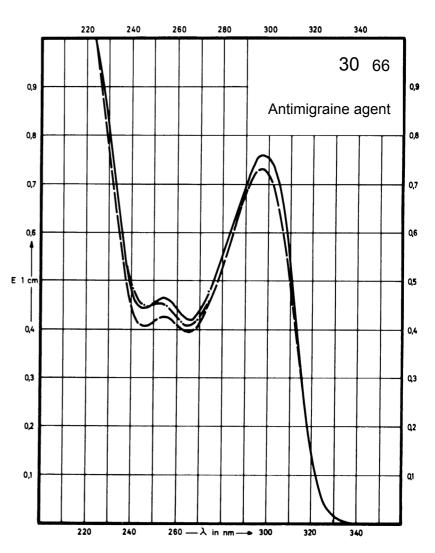

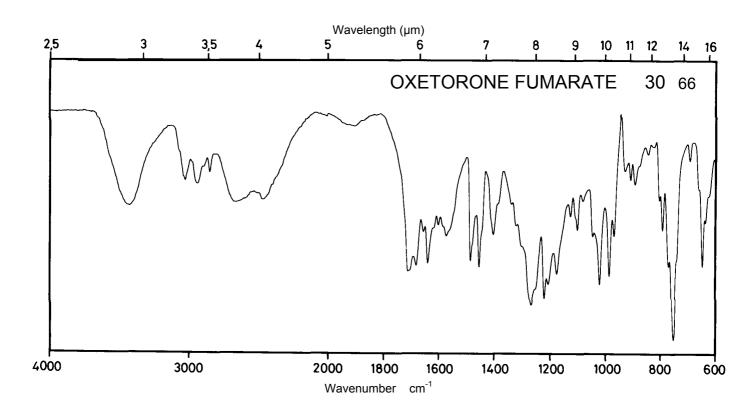

Name OBIDOXIME CHLORIDE


 M_{r}


359.2

Concentration

0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	380 nm 296 nm		285 nm	Decom- position observed
E ^{1%} 1cm	160 875		940	
ε	5750 31400		33700	

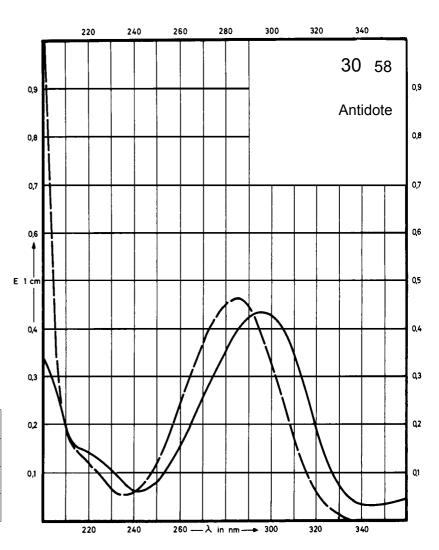

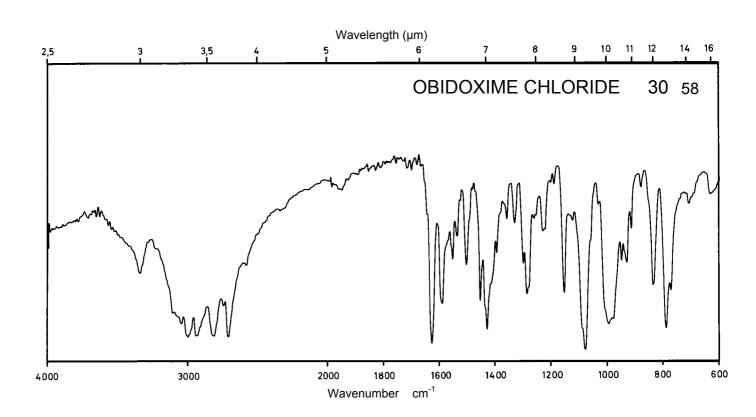







156.1


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm		282 nm	285 nm
E ^{1%} 1cm	370		370	293
3	5780		5780	4570

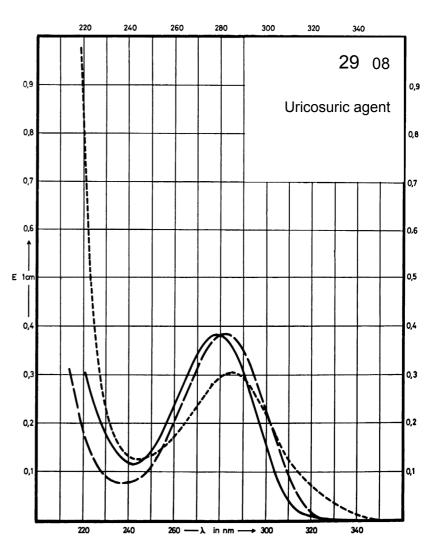

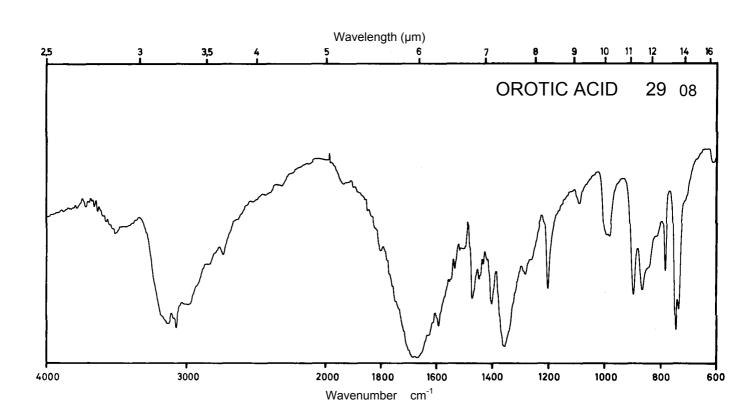




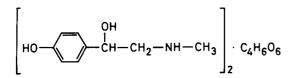


 M_r 484.5


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 225 nm		273 nm 222 nm	289 nm 241 nm
E ^{1%} 1cm	65 398		54 367	99 556
ε	3130 19300		2600 17790	4780 29900

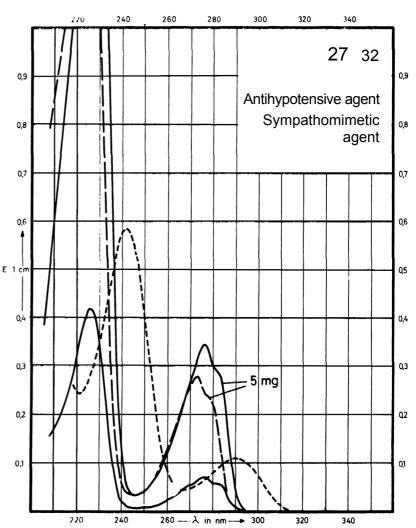

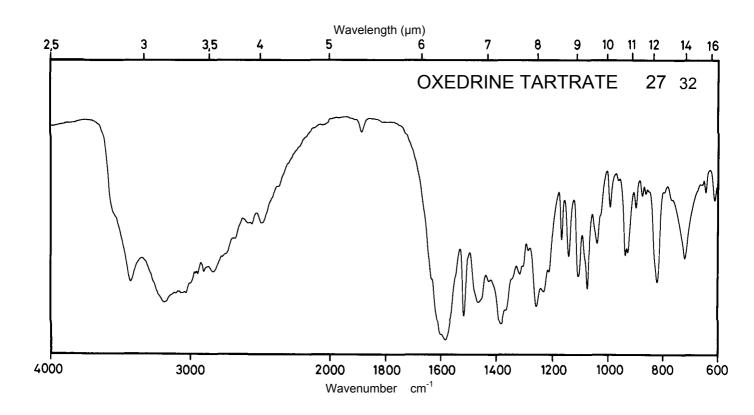


Name OXYMETAZOLINE HYDROCHLORIDE

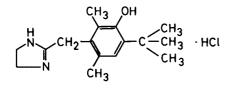


 M_r 296.8

Concentration 2 mg / 100 ml 6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm		280 nm	303 nm
E ^{1%} 1cm	80		60	144
ε	2310		1760	4260

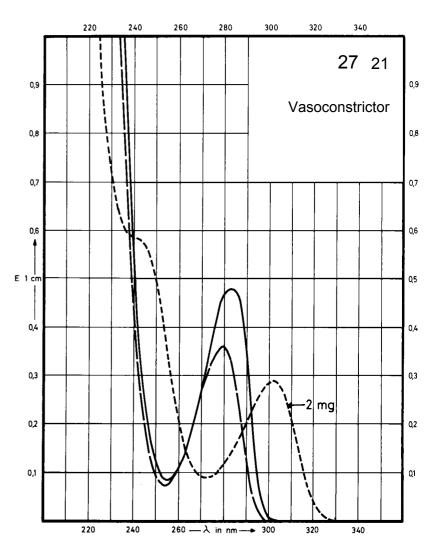

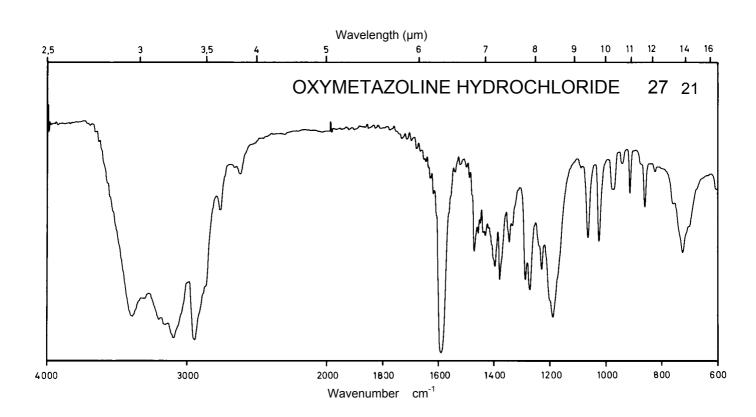

Name ORCIPRENALINE SULFATE


$$\begin{bmatrix} OH & CH_3 \\ -CH-CH_2-NH-CH \\ OH & CH_3 \end{bmatrix}_2 \\ \cdot H_2SO_4$$


 M_{r}

520.6


Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm		276 nm	296 nm
E ^{1%} 1cm	75		68	123
ε	3900		3530	6420

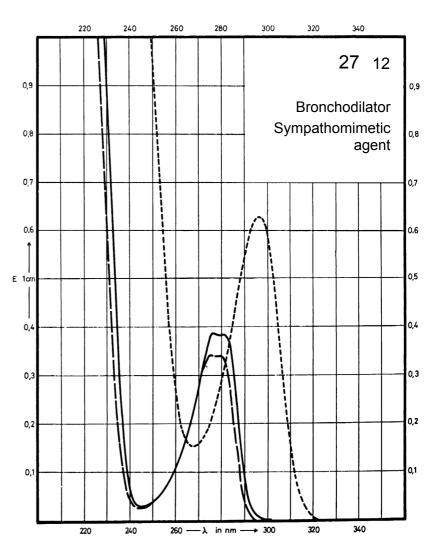

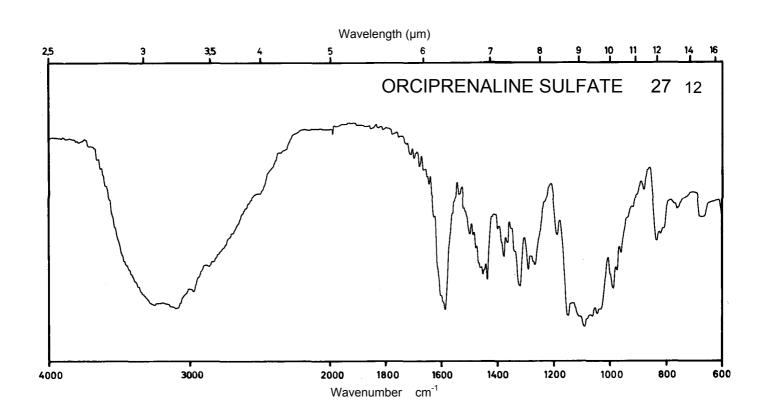




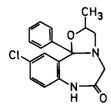


 M_r 328.8


Concentration 0.5 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm		368 nm 282 nm 239 nm	Decom- position observed
E ^{1%} 1cm	417		116 293 967	
3	13700		3830 9600 31800	

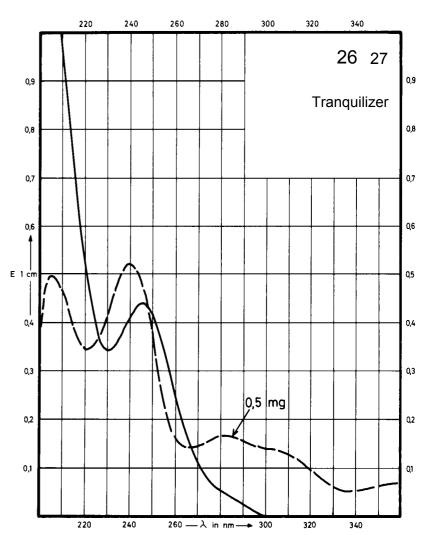

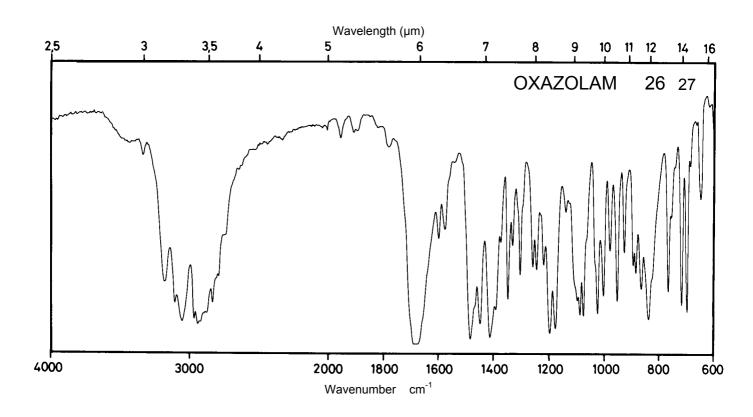




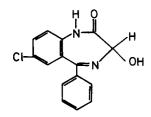


 M_r 286.7


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm		360 nm 283 nm	344 nm
E ^{1%} 1cm	84		120 400	102
3	2410		3440 11470	2920

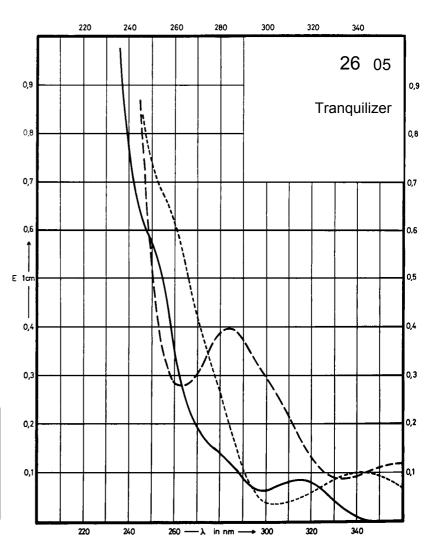

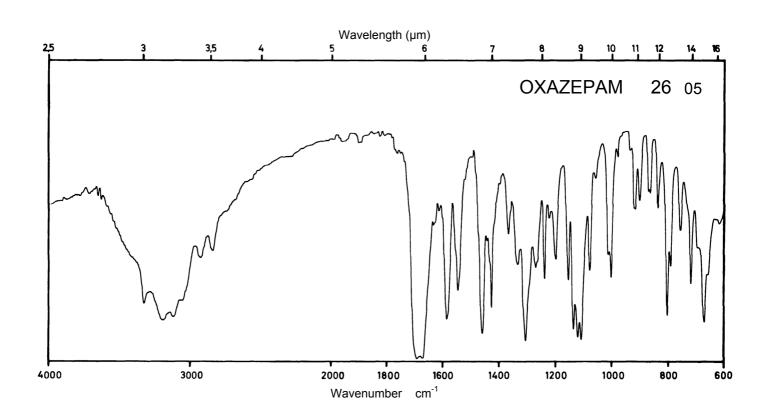




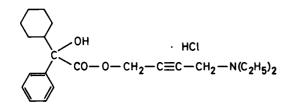


 M_r 394.0


Concentration 5 mg / 100 ml 100 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 264 nm 258 nm 252 nm 264 nm 258 nm 252 nm Maximum of absorption 4.18 5.38 4.68 5.80 E^{1%} 1cm 475 165 211 187 ε 228

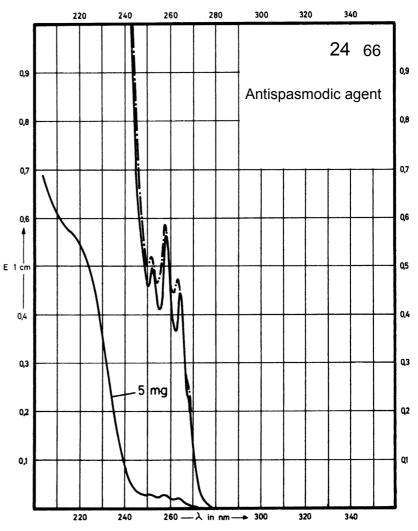

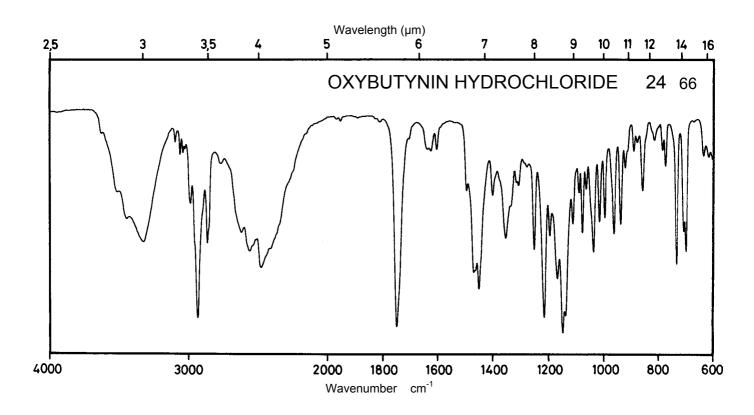



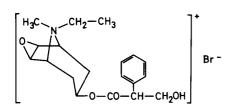



 M_r

412.3

Concentration 150 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm		257 nm	
E ^{1%} 1cm	4.60		4.68	
ε	190		195	

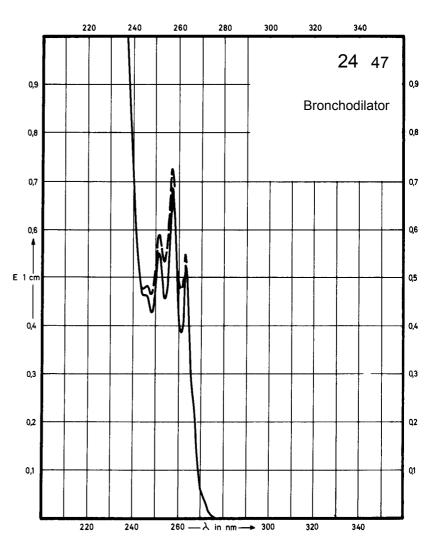

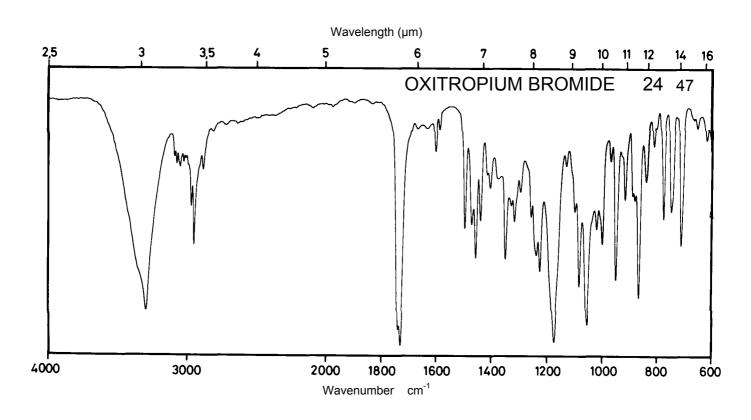





 M_r

461.5

Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm		264 nm	
E ^{1%} 1cm	12.0		13.7	
ε	550		630	

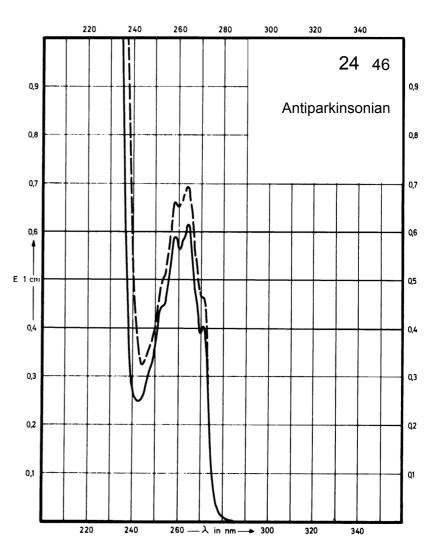

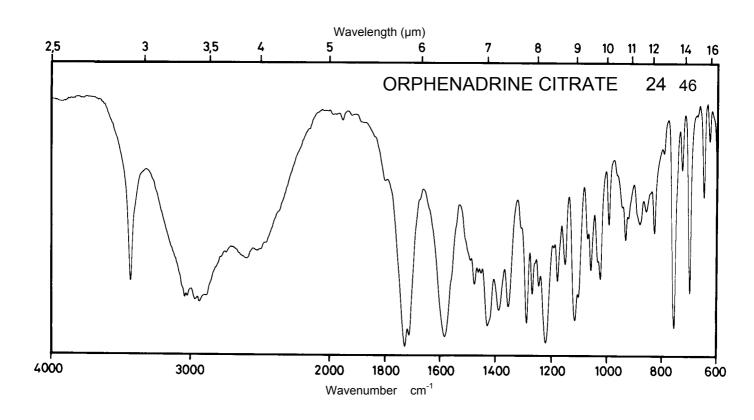




 $C_{43}H_{66}N_{12}O_{12}S_2$


 M_{r}

1007.2

Concentration 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm		274 nm	293 nm
E ^{1%} 1cm	13.4		11.9	19.0
ε	1350		1200	1910

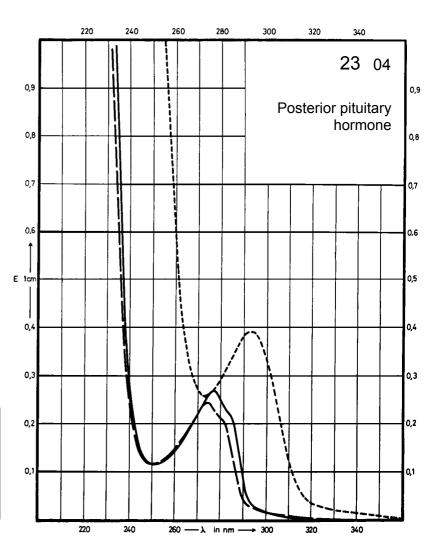

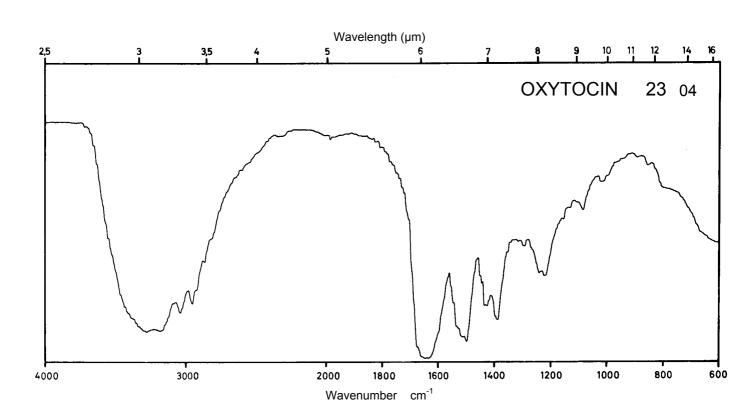



Name OXABOLONE CIPIONATE

 M_{r} 414.6


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm			
E ^{1%} 1cm	318			
3	13180			

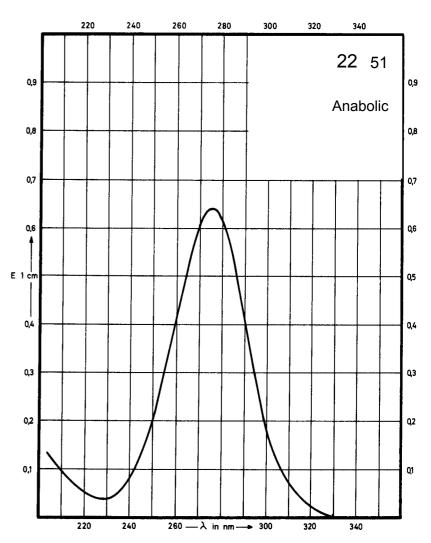

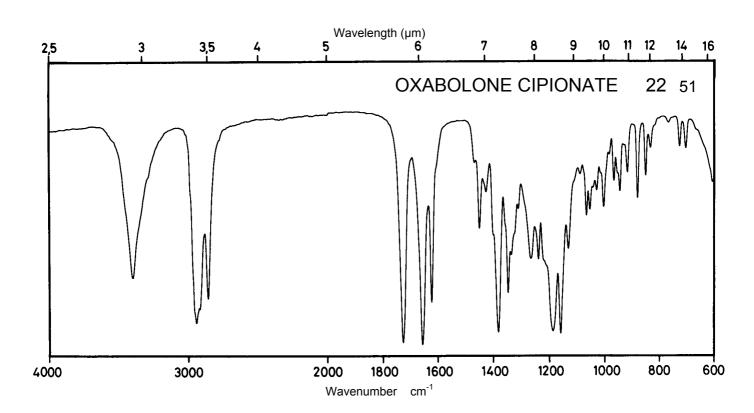

Name OXYMETHOLONE

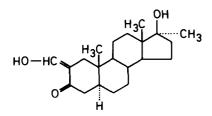



 M_{r}


332.5


Concentration 1 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm			315 nm
E ^{1%} 1cm	244			543
3	8130			18060

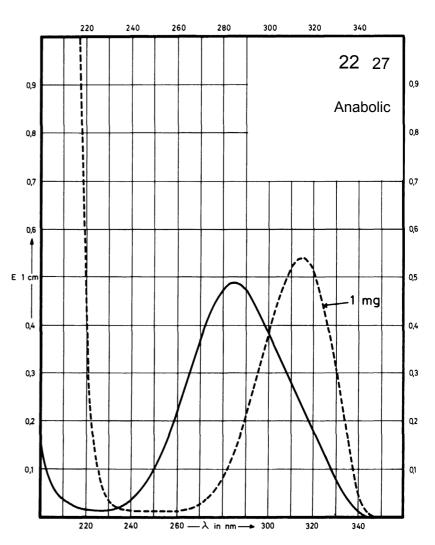


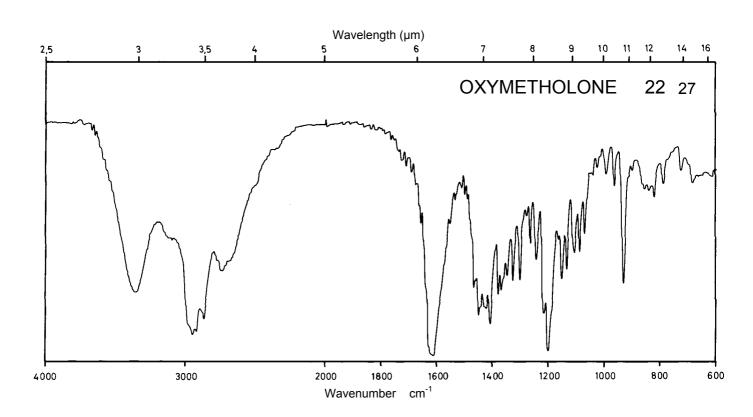

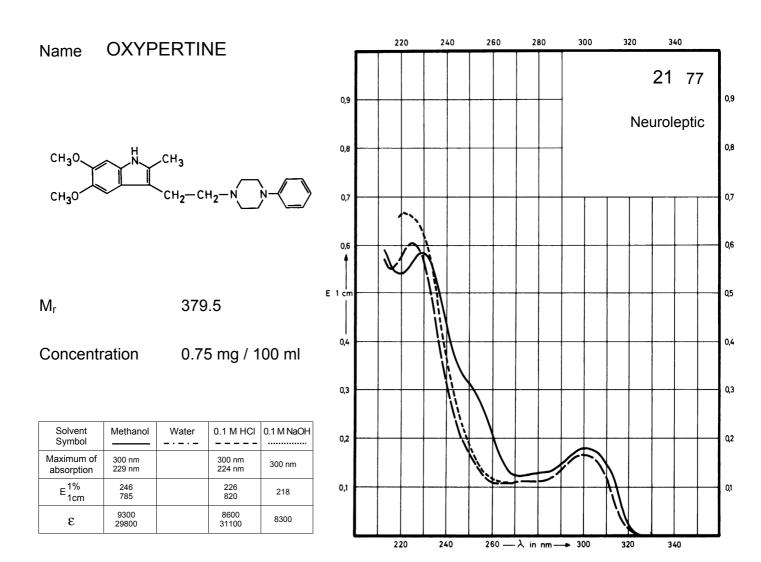

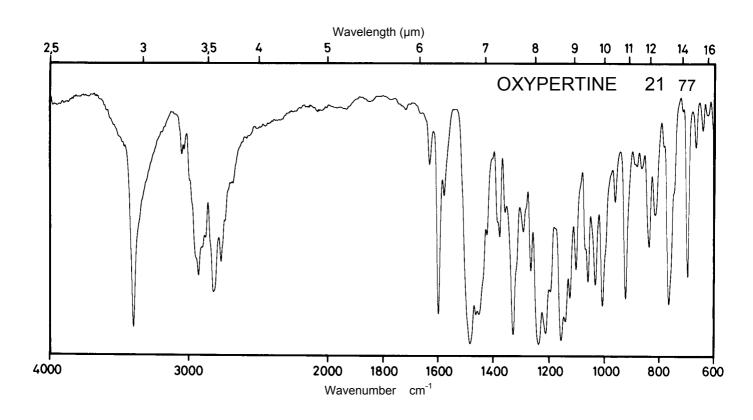


Name OPIPRAMOL DIHYDROCHLORIDE

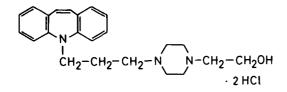


 M_{r}


436.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	256 nm		253 nm	253 nm
E ^{1%} 1cm	760		665	665
ε	33170		29020	29020

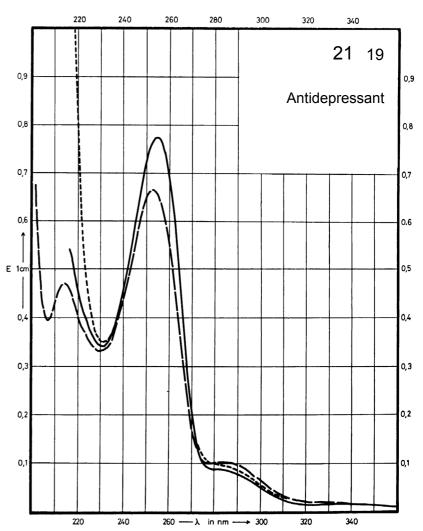

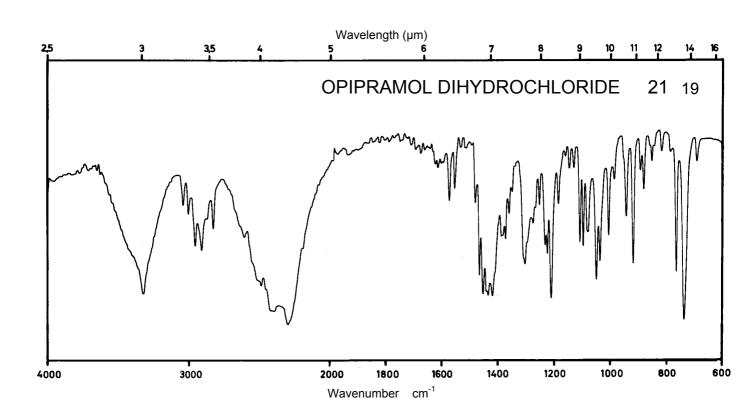



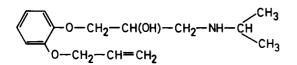



301.8


Concentration

2.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		272 nm	272 nm
E ^{1%} 1cm	77		77	80
ε	2330		2330	2410

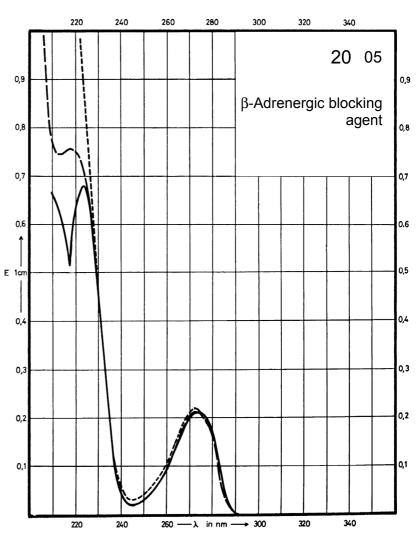


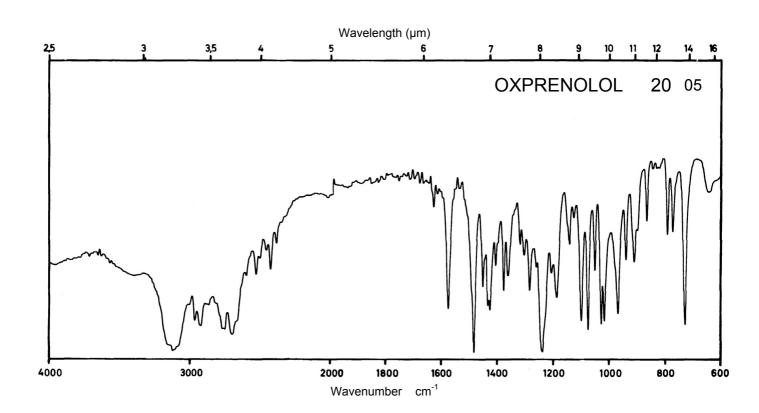


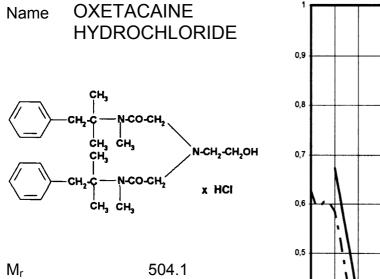


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm			
E ^{1%} 1cm	7.2 10.6 11.0			
ε	360 530 550			

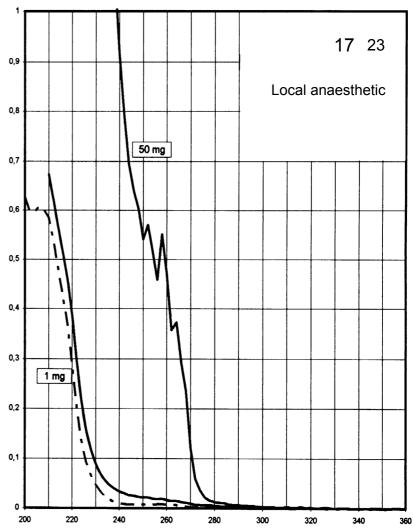

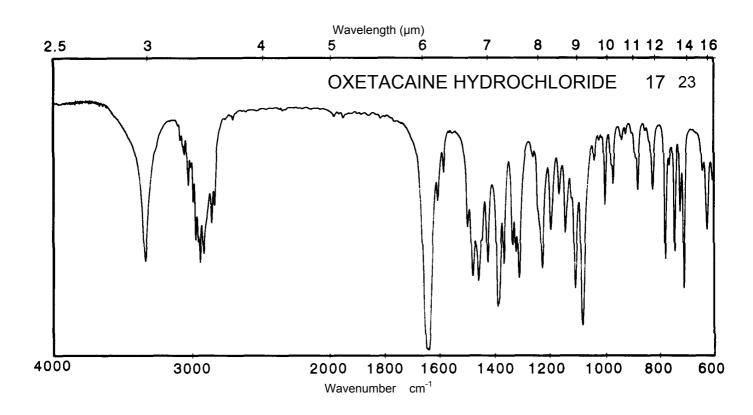


Name OXYBUPROCAINE HYDROCHLORIDE

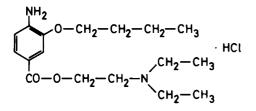



 M_r


344.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm 232 nm		297 nm 235 nm	306 nm
E ^{1%} 1cm	520 375		107 238	Decom- position observed
ε	17900 12900		3700 8200	

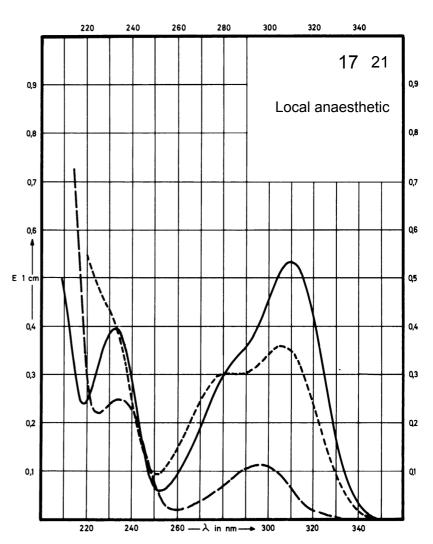


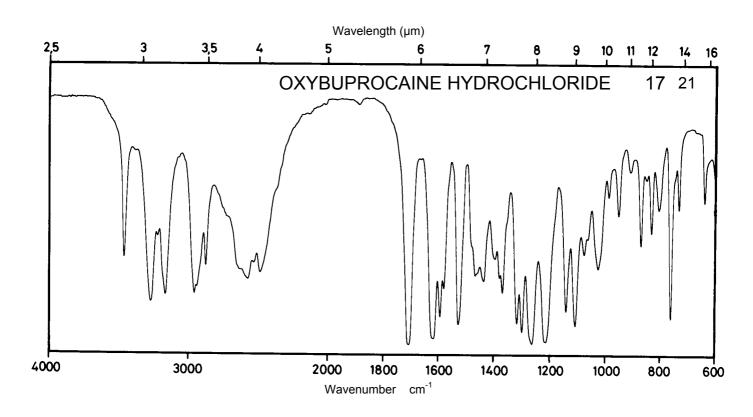

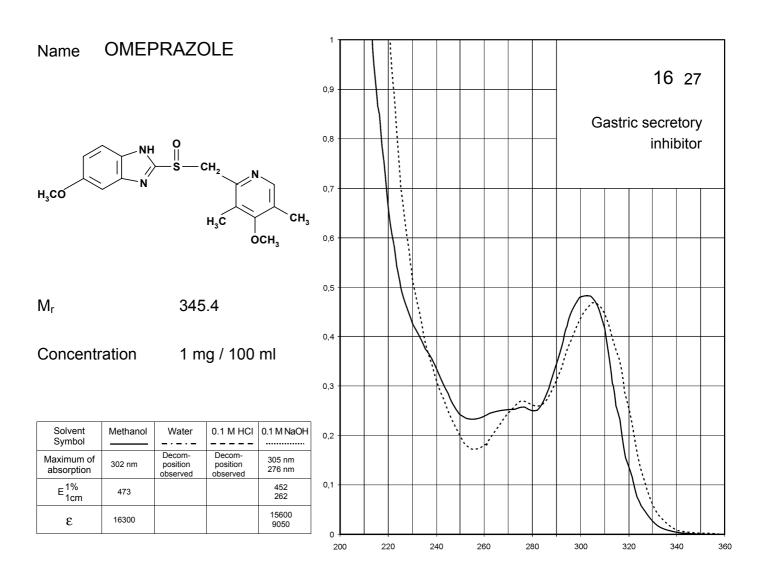

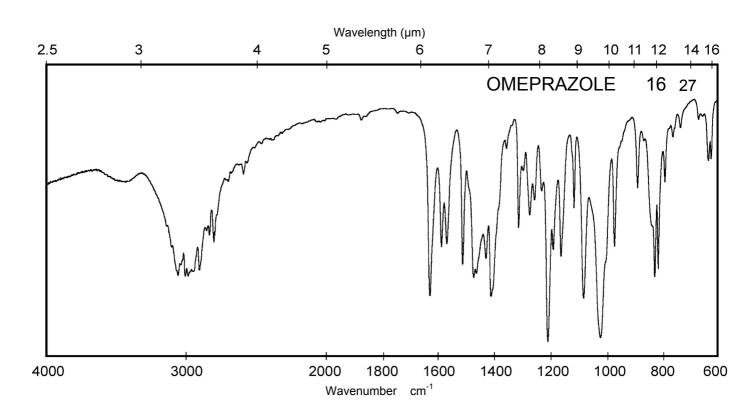


Name OCTOPAMINE HYDROCHLORIDE

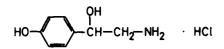


 M_{r} 189.6


Concentration 1.25 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 277 nm 225 nm 290 nm 241 nm 273 nm 273 nm absorption E^{1%} 1cm 85 500 70 70 1600 9500 2400 13400 ε 1300 1300

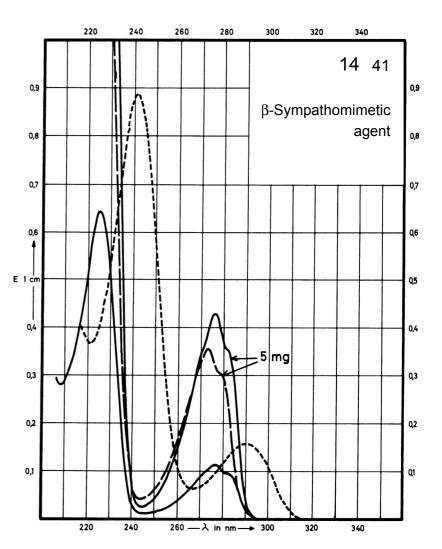

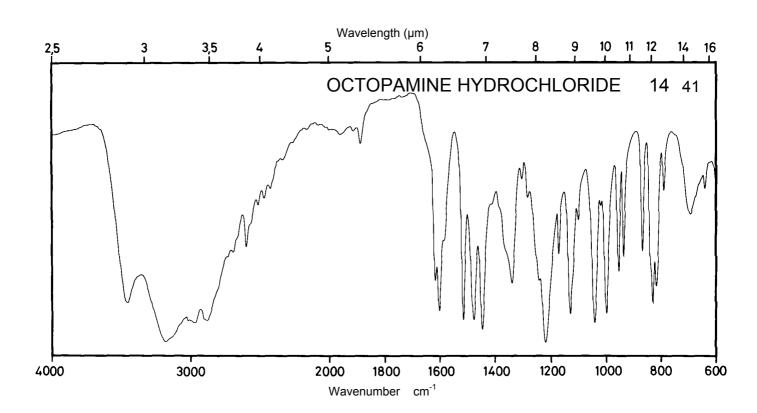




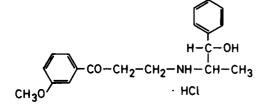


 M_r

349.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	309 nm 252 nm		310 nm 253 nm	
E ^{1%} 1cm	74 266		69 274	
ε	2580 9310		2420 9570	

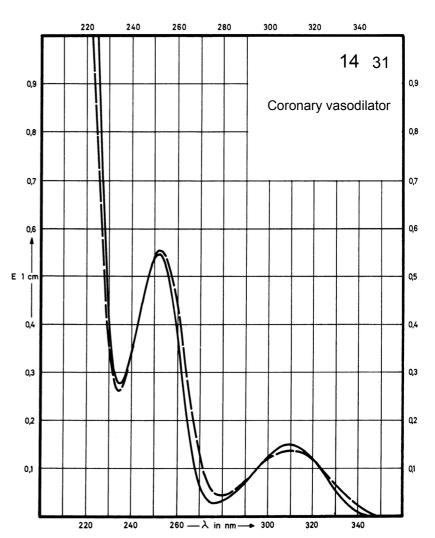

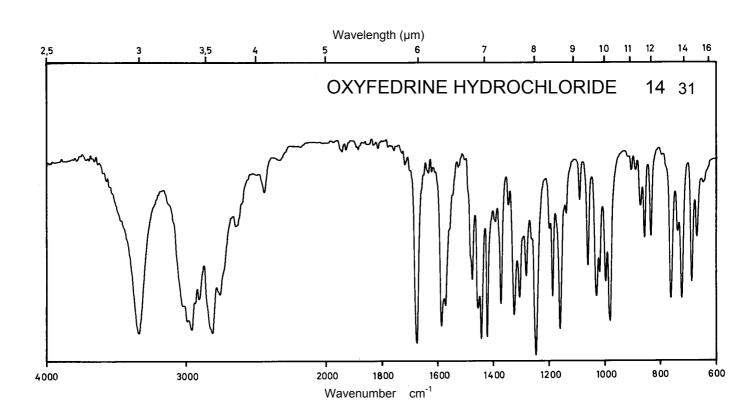







145.2


Concentration 0.25 mg / 100 ml 2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Gyllibol				
Maximum of	310 nm		358 nm	344 nm
	241 nm		318 nm	334 nm
absorption	24 1 11111		251 nm	252 nm
10/-	178		118	190
E ^{1%} 1cm	2540		112	200
1cm	2540		3090	2170
	0500		1710	2760
3	2580		1630	2900
C	36880		44870	31510

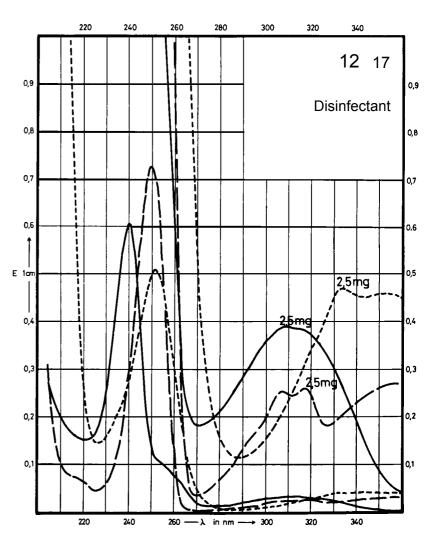

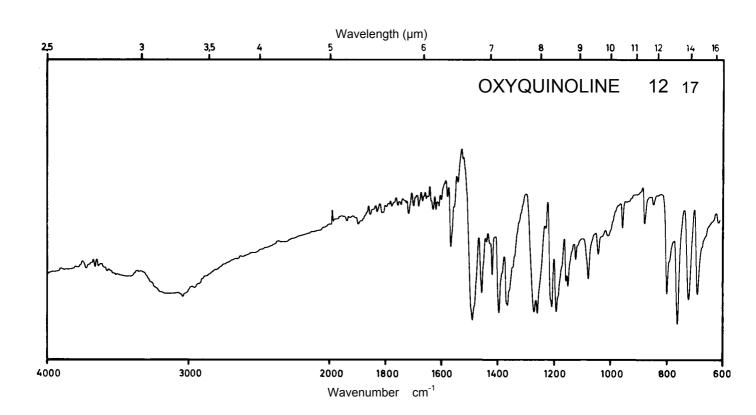



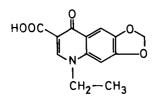



261.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			336 nm 322 nm 259 nm	340 nm 266 nm 258 nm
E ^{1%} _{1cm}			328 333 1660	457 1320 1250
ε			8600 8700 43400	11900 34500 32700

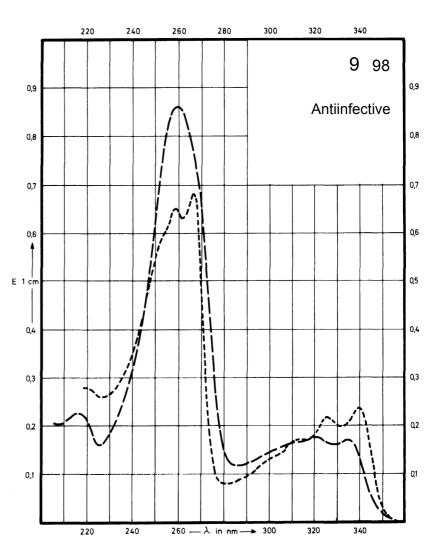

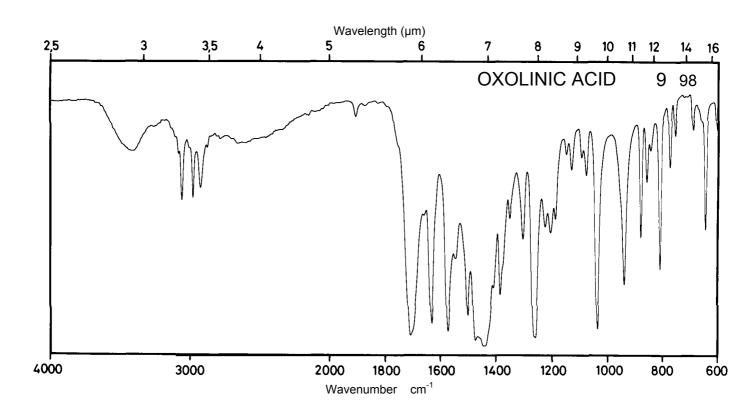



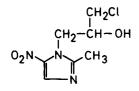



219.6

Concentration

1.2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm		277 nm	320 nm
E ^{1%} 1cm	392		281	400
ε	8620		6180	8780

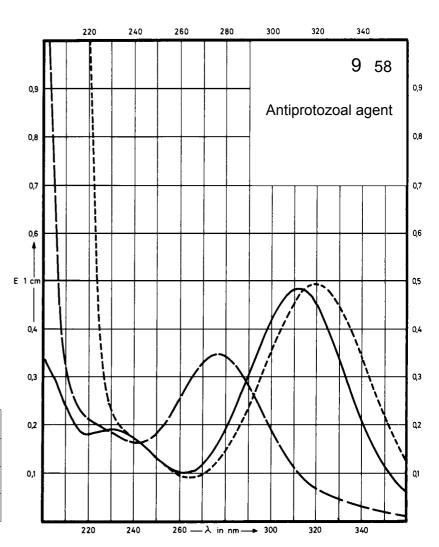

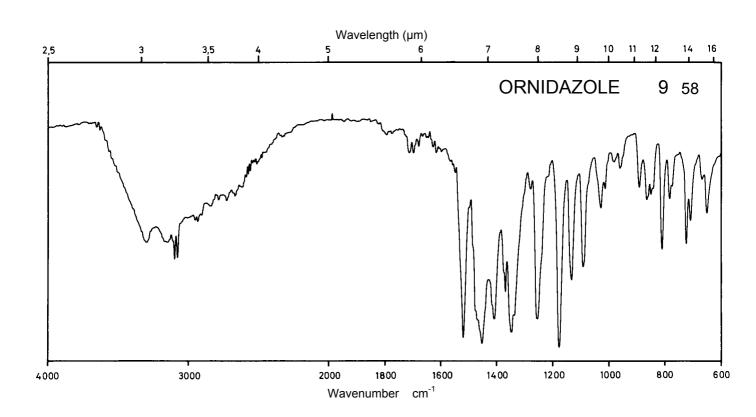


Name OXOLAMINE CITRATE


 M_r

Concentration 2 mg / 100 ml 40 mg / 100 ml


437.5

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 285 nm 277 nm 238 nm 11.7 16.2 Maximum of 239 nm 239 nm 235 nm absorption 270 272 300 300 510 710 13100 ε 11800 11900 13100

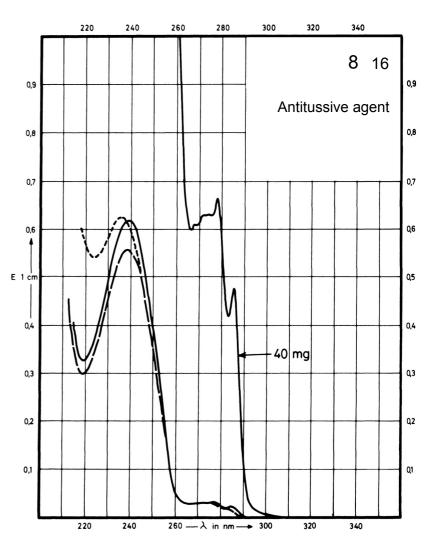

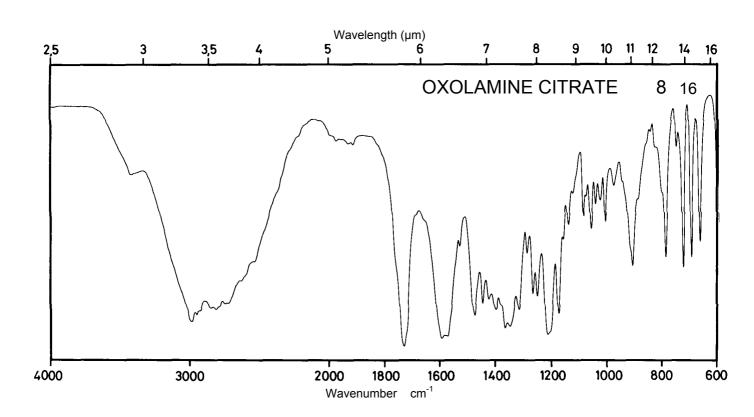




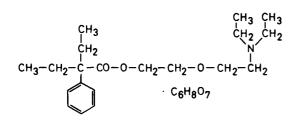



 M_r


527.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
-,	005		005	
Maximum of	265 nm		265 nm	
	259 nm		259 nm	
absorption	253 nm		253 nm	
1%	2.62		2.64	
E ^{1%} 1cm	3.49		3.55	
1cm	2.83		2.91	
	138		139	
3	184		187	
C	1/10		154	

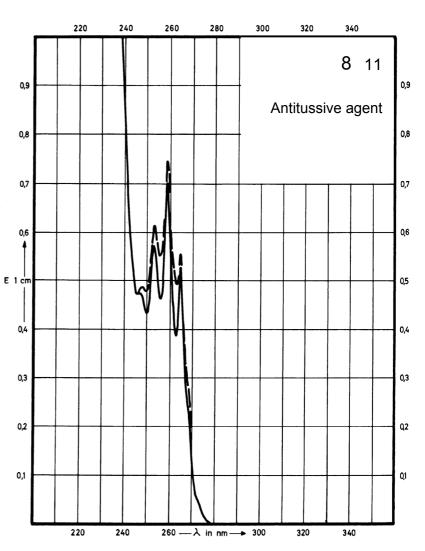


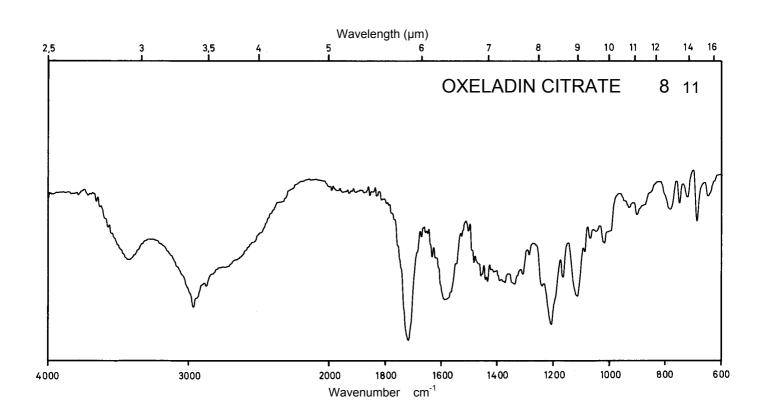

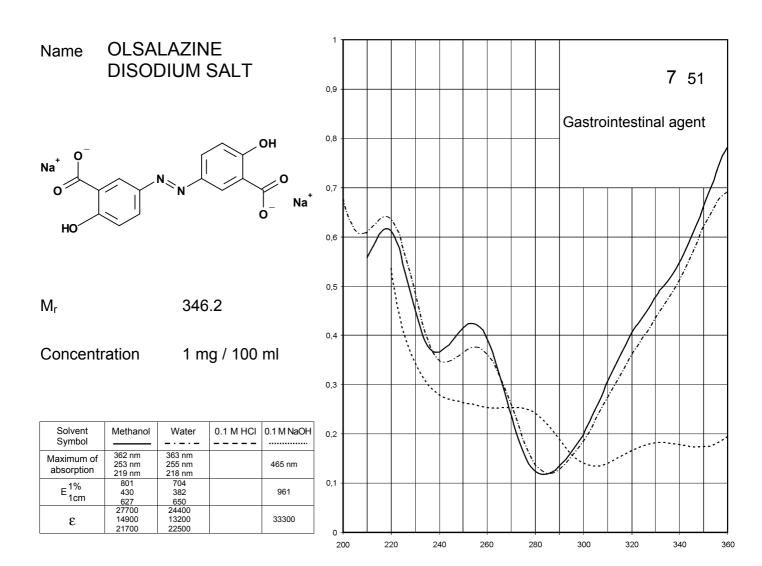

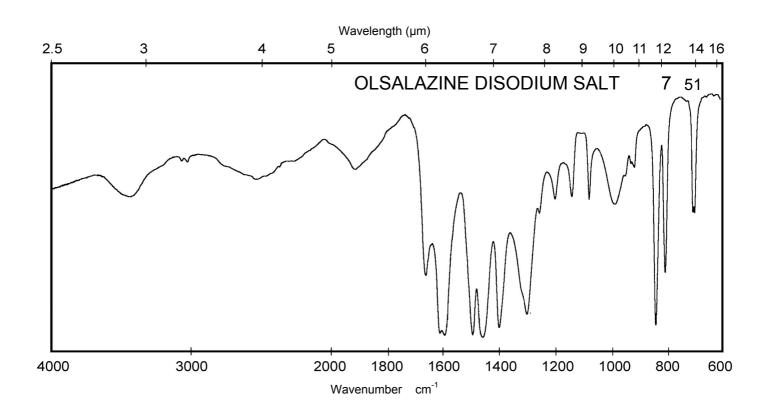





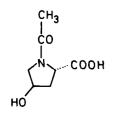


173.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

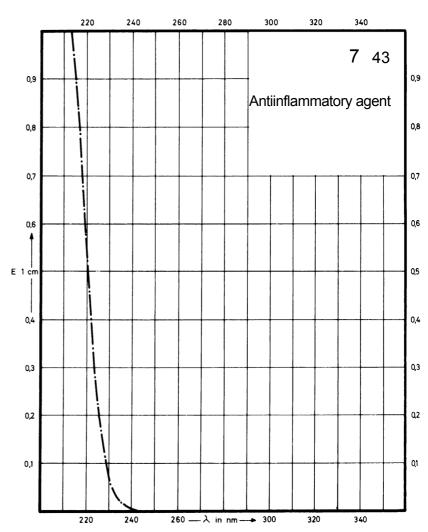

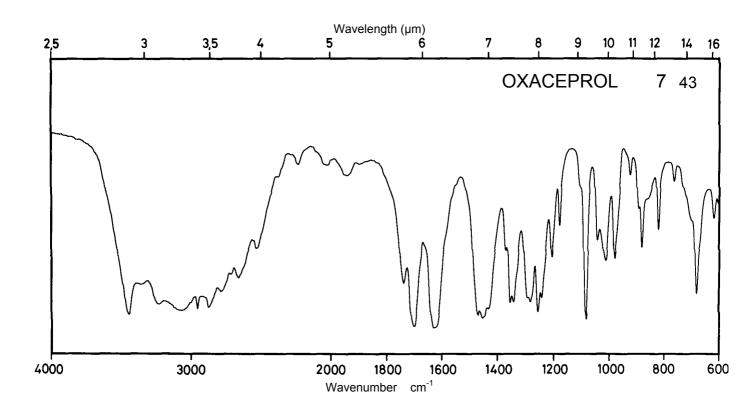





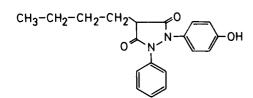


324.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm		236 nm	254 nm
E ^{1%} 1cm	480		476	702
3	15570		15440	22770

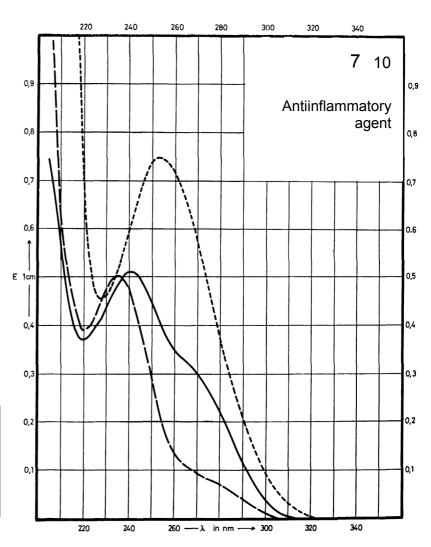

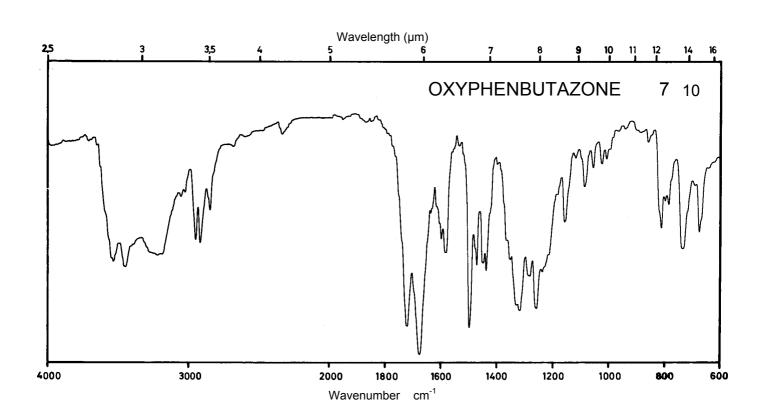



Name OXYTETRACYCLINE HYDROCHLORIDE

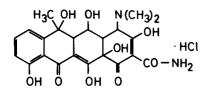


 M_r 496.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	360 nm 266 nm		353 nm 267 nm	376 nm 269 nm
E 1% 1cm	293 364		276 380	331 320
ε	14540 18110		13730 18880	16460 15900

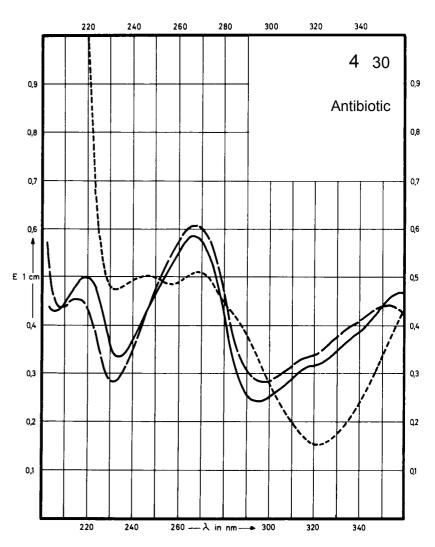

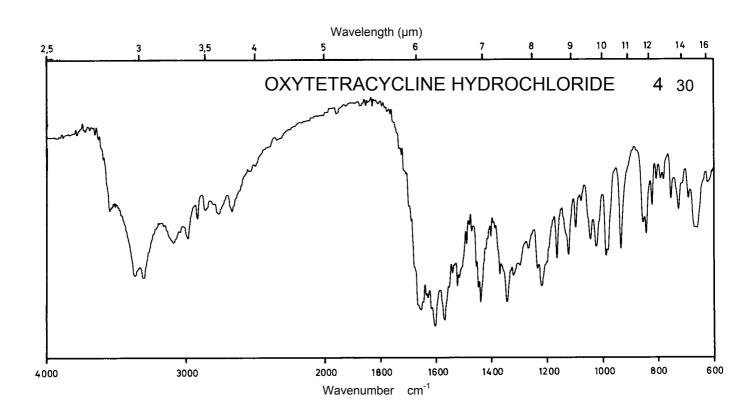





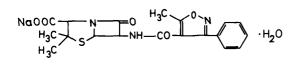

441.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

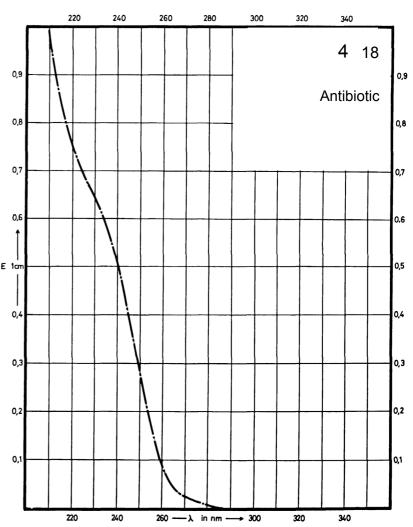

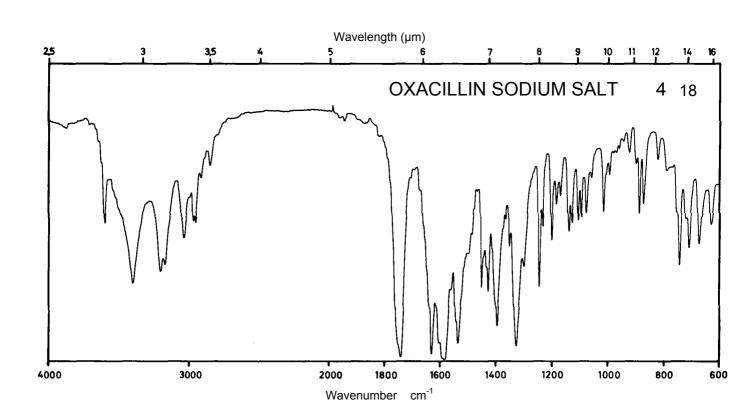



Name OXYTETRACYCLINE

 M_{r}


496.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm		352 nm 266 nm	269 nm
E ^{1%} 1cm	335		275 378	327
3	16630		13650 18770	16240

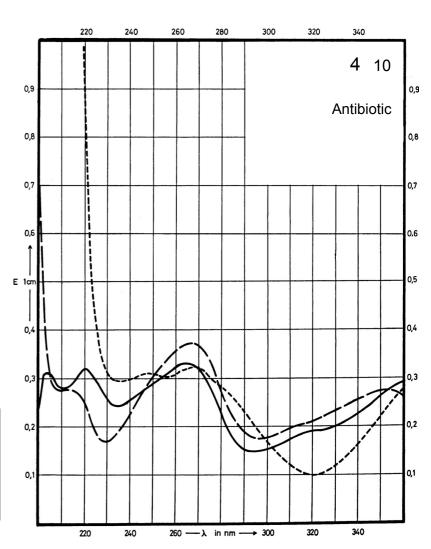

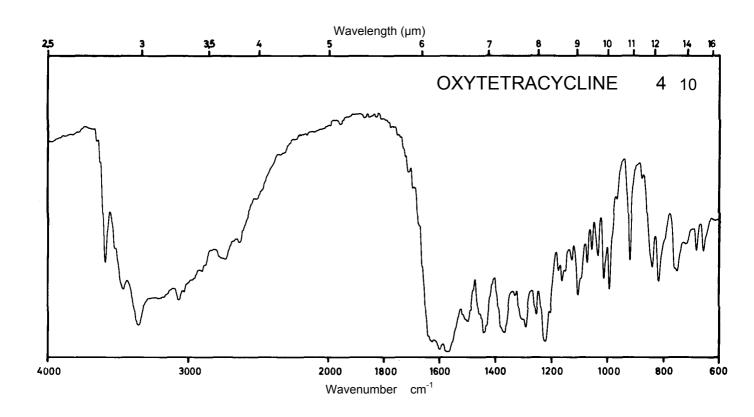


Name ONDANSETRON HYDROCHLORIDE

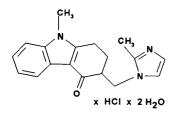


 M_{r}


365.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	303 nm	310 nm	310 nm	310 nm
	267 nm	267 nm	267 nm	267 nm
	246 nm	249 nm	249 nm	248 nm
E ^{1%} 1cm	373 338 462	397 360 445	404 360 442	400 362 447
ε	13600	14500	14800	14600
	12400	13200	13200	13200
	16900	16300	16200	16300

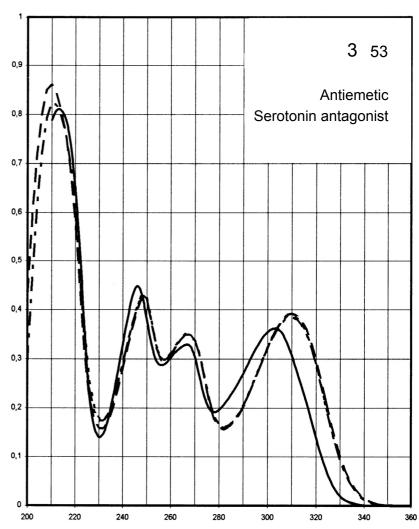

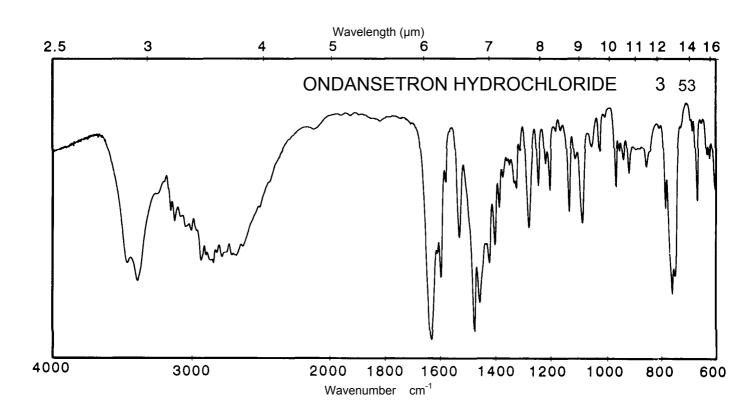





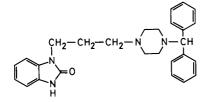


 M_{r} 426.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm 225 nm		278 nm 225 nm	287 nm
E ^{1%} 1cm	158 386		153 420	170
ε	6700 16500		6500 17900	7200

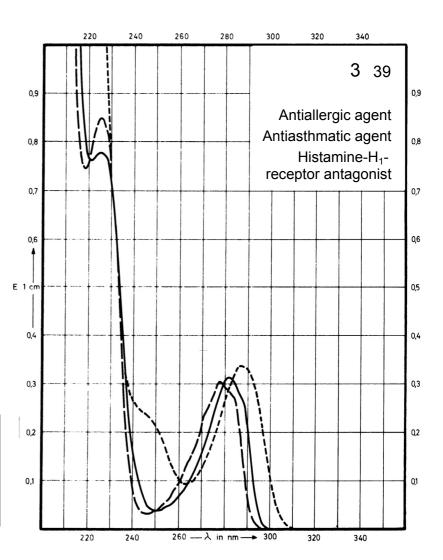

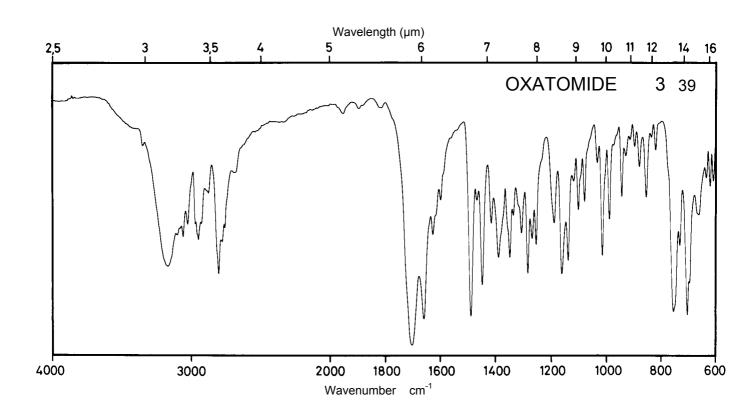


Name OXYPENDYL DIHYDROCHLORIDE

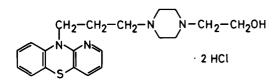


 M_r 443.4


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	320 nm 252 nm		372 nm 242 nm	320 nm 251 nm
E ^{1%} 1cm	111 668		113 547	106 581
ε	4940 29600		5000 24250	4700 25780

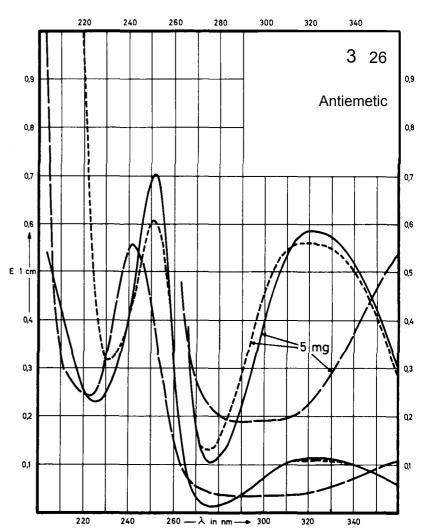

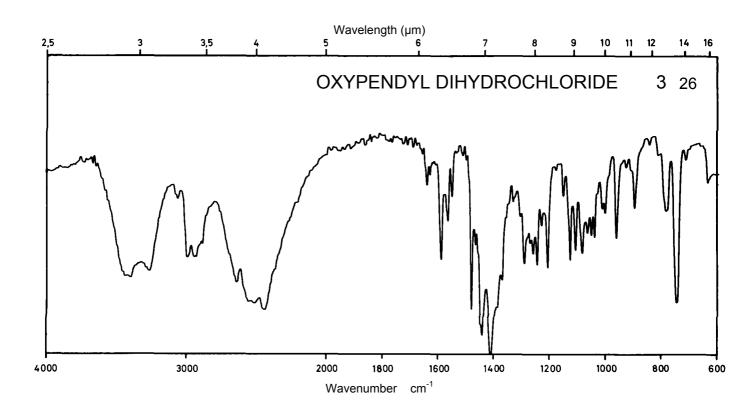

Name OXOMEMAZINE HYDROCHLORIDE

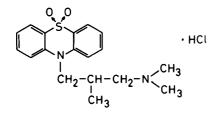


 M_{r}


366.9

Concentration


0.75 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	332 nm 268 nm 226 nm		333 nm 267 nm 226 nm	335 nm 298 nm 270 nm
E ^{1%} 1cm	150 360 815		153 346 825	160 230 353
ε	5500 13210 29900		5610 12690 30270	5870 8440 12950

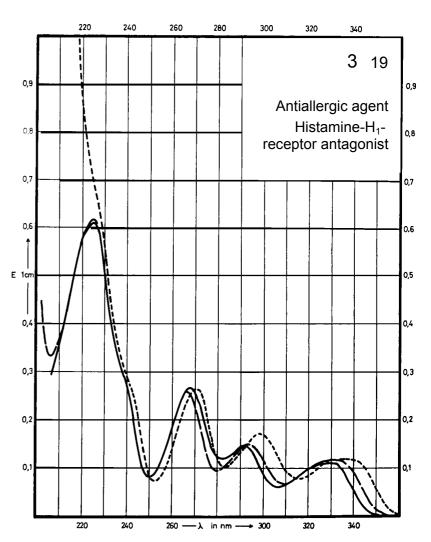

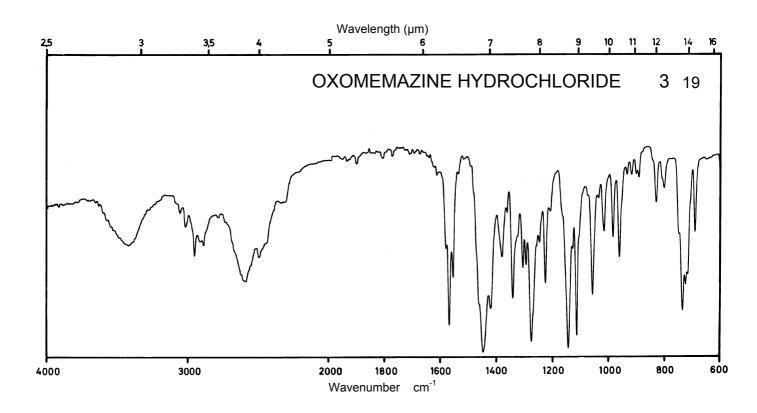





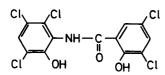


401.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	301 nm			324 nm
E ^{1%} 1cm	209			338
ε	8390			13570

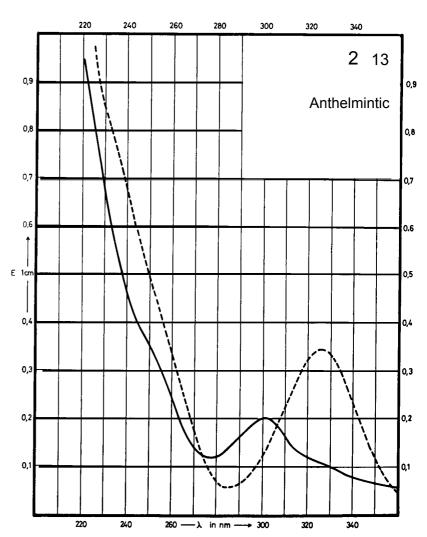

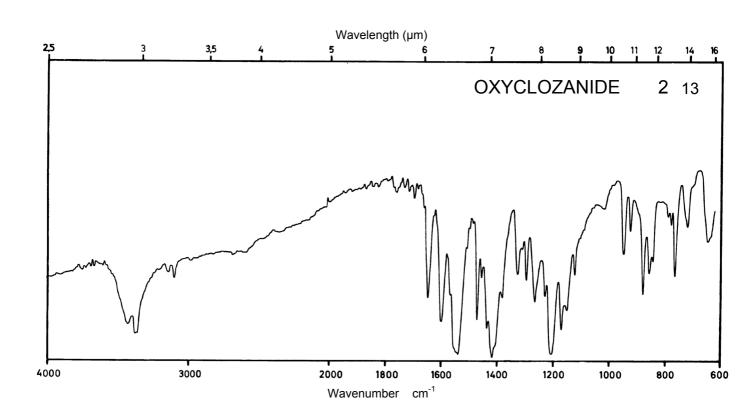





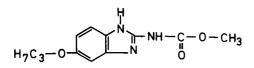


249.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm		290 nm	305 nm 254 nm
E ^{1%} 1cm	544		590	767 462
ε	13560		14710	19120 11520

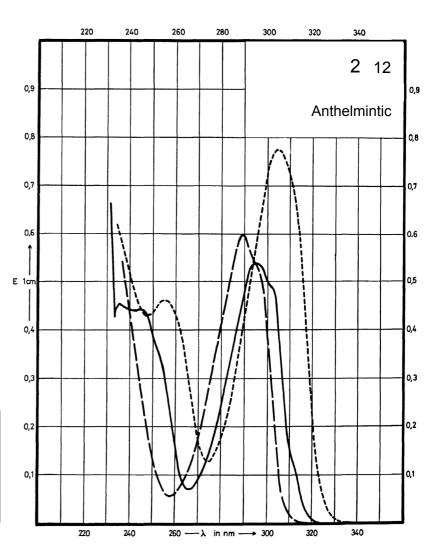

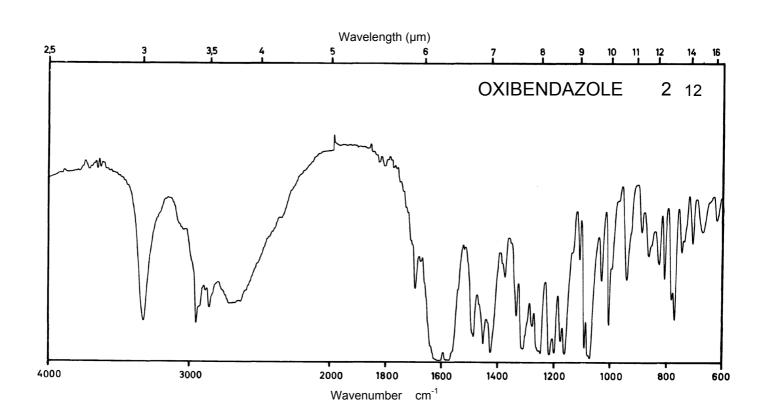



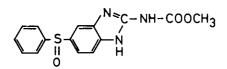



315.3


Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm		290 nm 283 nm	
E ^{1%} 1cm	570		600 550	
ε	17970		18920 17340	

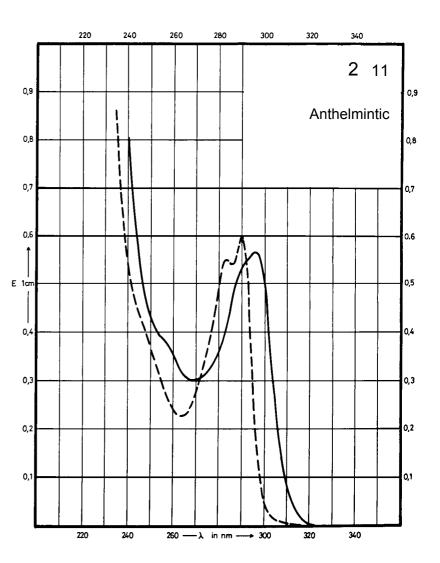


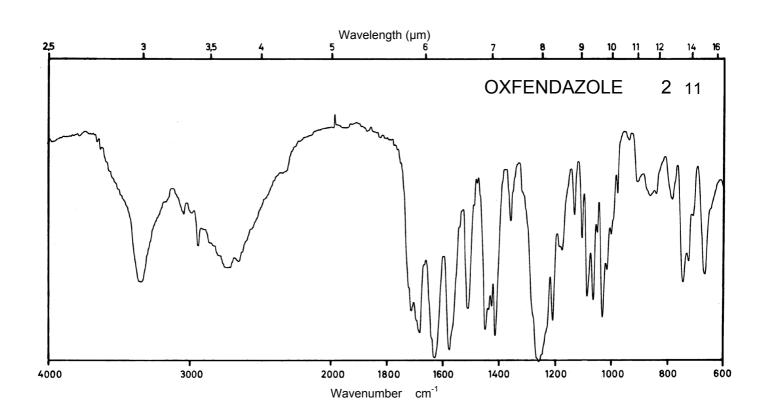

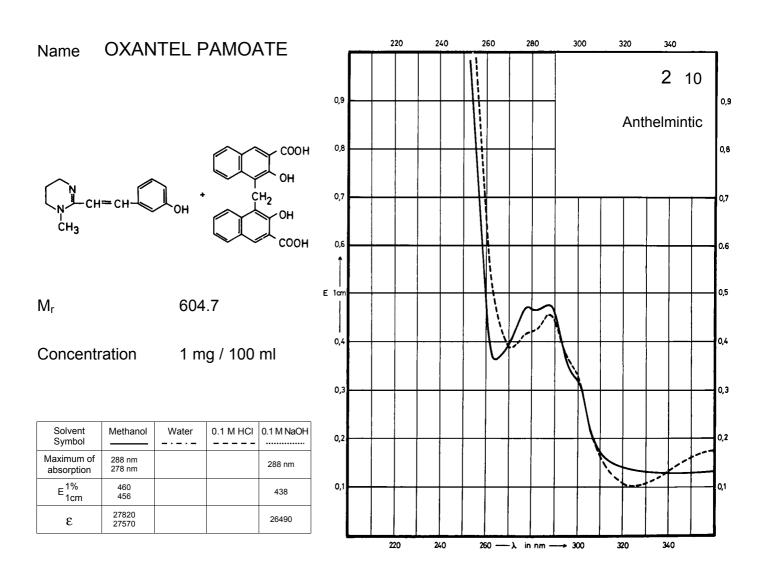

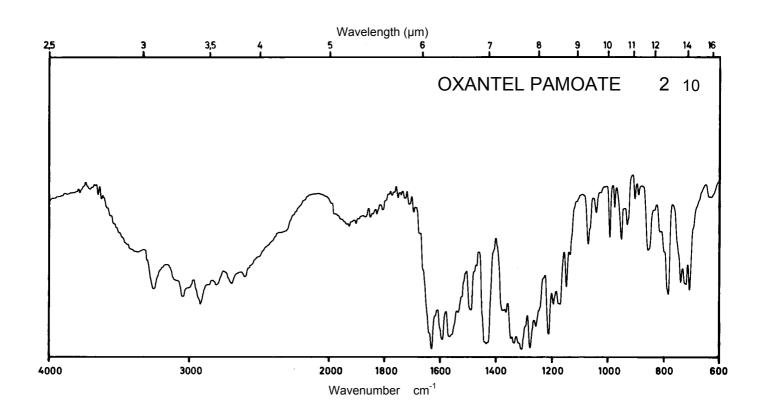





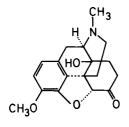


 M_r 315.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm			
E ^{1%} 1cm	37.6			
ε	1180			

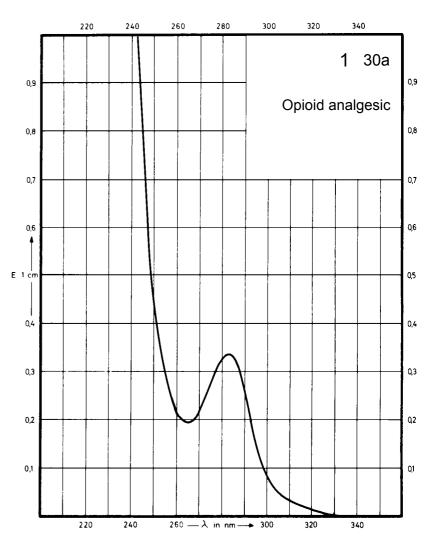

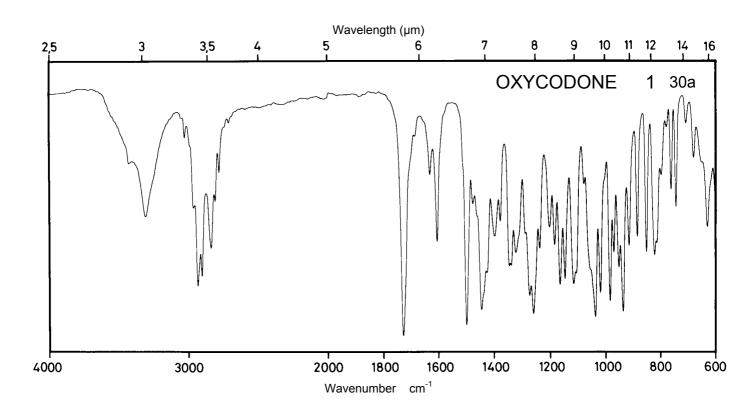


Name OXYCODONE HYDROCHLORIDE

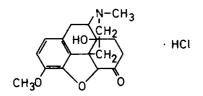


 M_{r}


351.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm		280 nm	280 nm
E ^{1%} 1cm	38.6		35.3	33.7
ε	1360		1240	1190

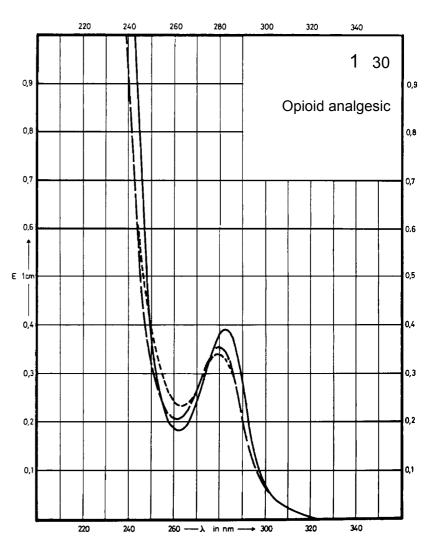

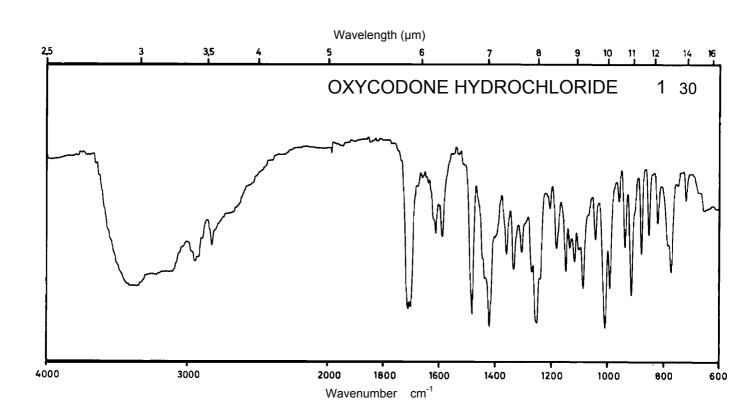




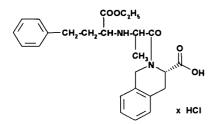


 M_r

475.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm			
E ^{1%} _{1cm}	9.5 9.6			
3	450 450			

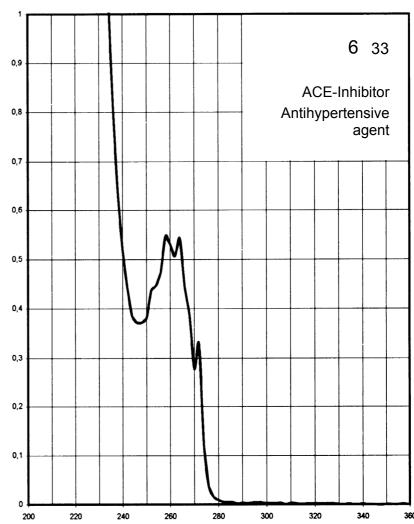

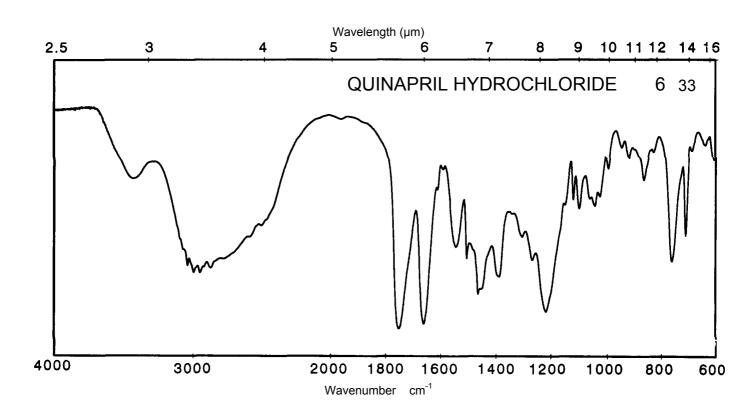



Name QUINIDINE

 M_{r}


324.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	332 nm 230 nm		346 nm 249 nm	330 nm
E ^{1%} 1cm	163 1022		176 903	146
ε	5290 33150		5710 29290	4740

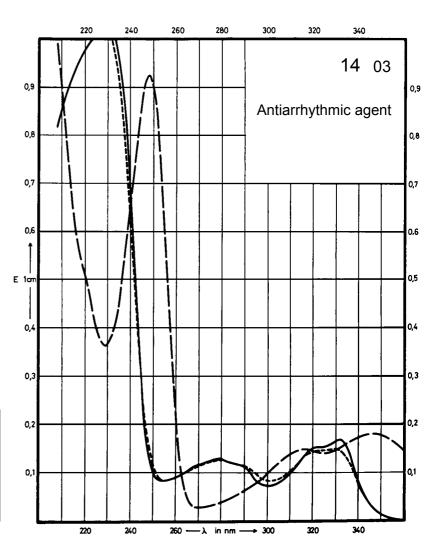

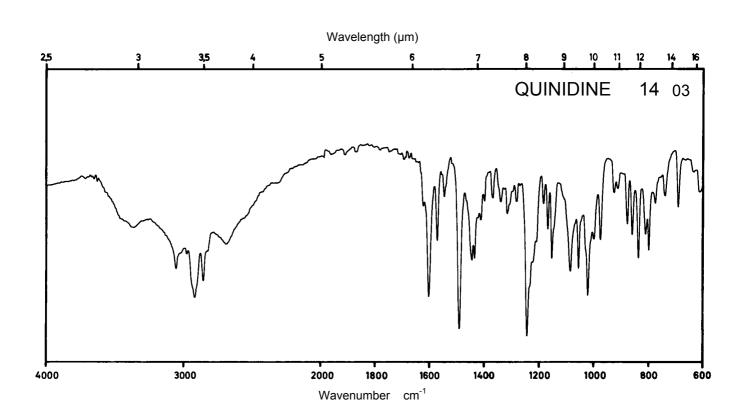




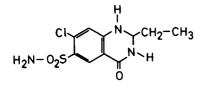


 M_{r} 289.7


Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 345 nm 280 nm 234 nm 77 328 345 nm 278 nm 234 nm 81 346 348 nm 268 nm 233 nm 72 238 Maximum of absorption E^{1%} 1cm 1525 2350 1470 2230 ε 10020 9500 6890 44180 42590

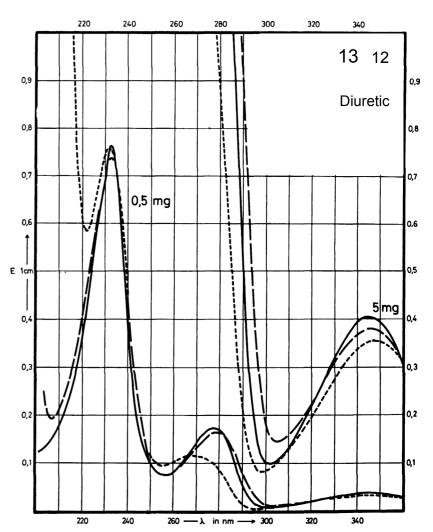

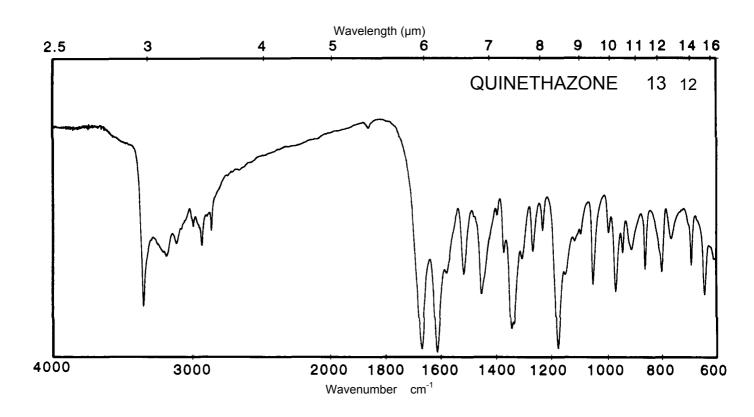




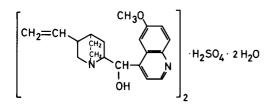



 $\,M_r\,$


782.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	333 nm 234 nm		346 nm 250 nm	330 nm 280 nm
E ^{1%} 1cm	130 797		137 737	117 92
ε	10100 62400		10700 57700	9200 7200

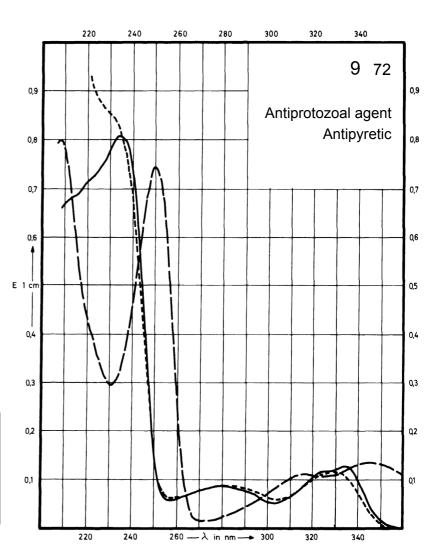

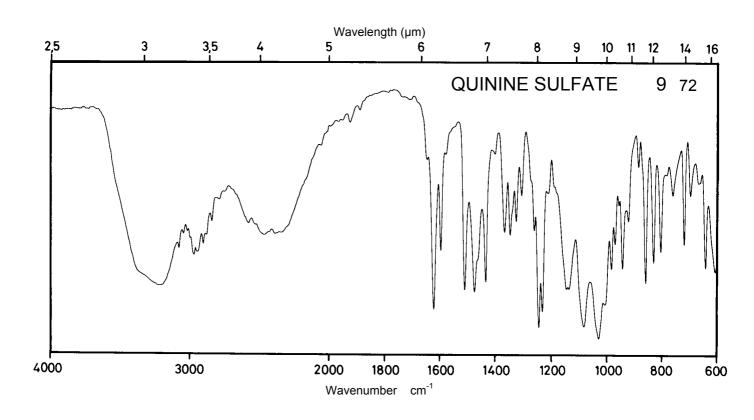

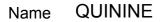

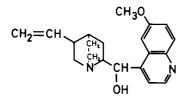



324.4

Concentration


0.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	333 nm 231 nm		346 nm 250 nm	330 nm 280 nm
E ^{1%} 1cm	156 1000		167 898	144 110
ε	5060 32400		5400 29100	4660 3570

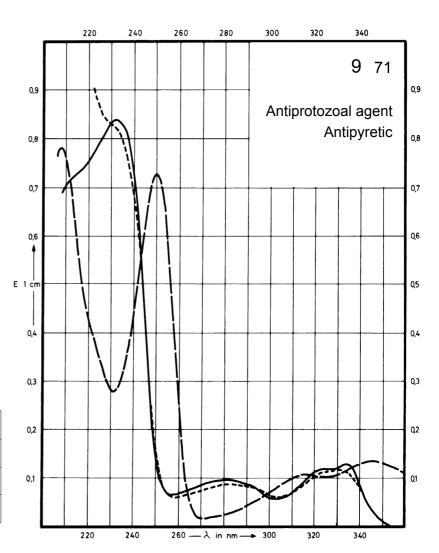

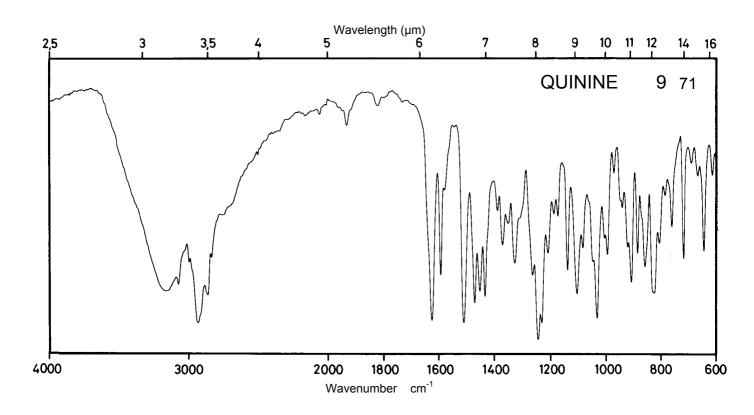





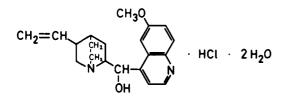



396.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	333 nm 234 nm		346 nm 250 nm	330 nm 280 nm
E ^{1%} 1cm	136 816		140 756	120 94
ε	5390 32380		5580 30000	4780 3730

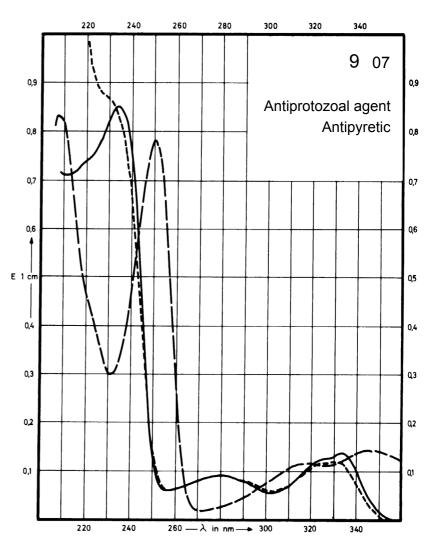


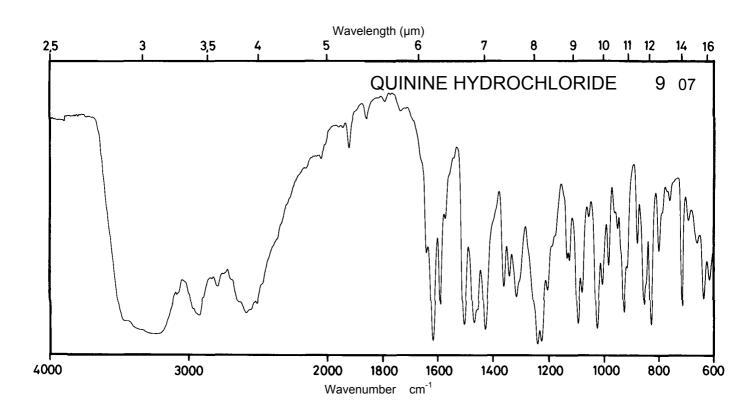

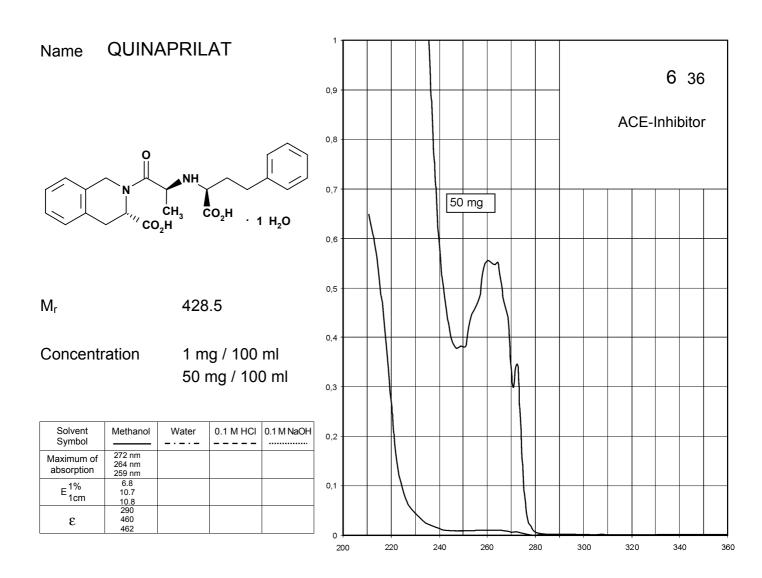

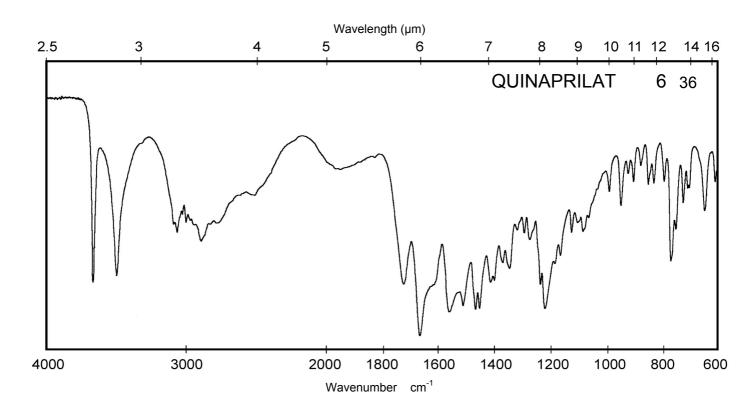





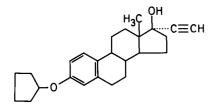


 M_r 364.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 280 nm			
E ^{1%} 1cm	48 53			
ε	1750 1950			

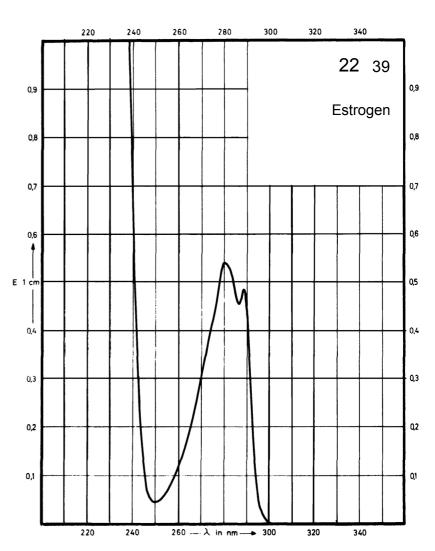

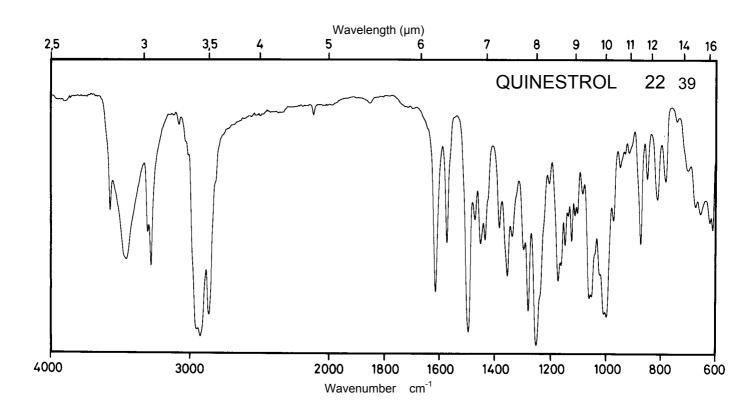




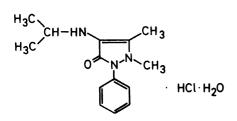


 M_r 299.8


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm 243 nm	260 nm	260 nm	271 nm
E ^{1%} *	332 270	318	328	291
ε	9360 7610	8960	9230	8211

^{*} Calculated on dried substance

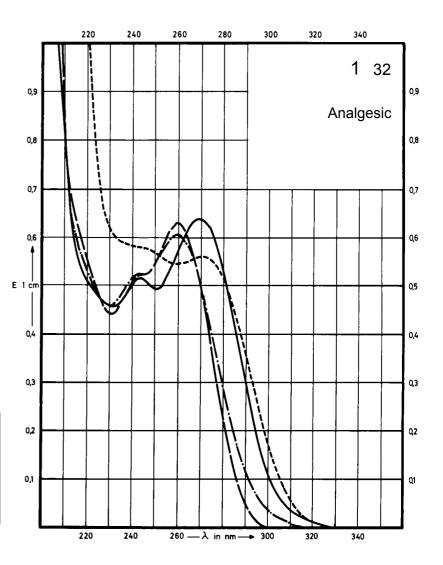

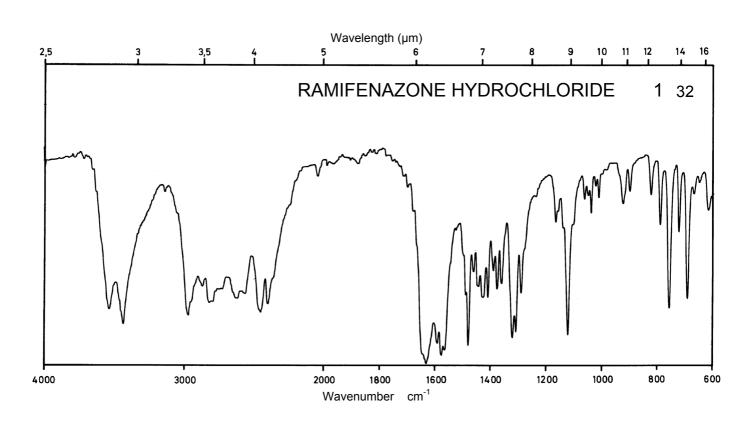




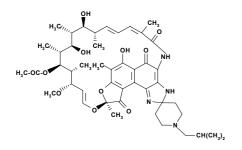



 M_{r}


847.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm	319 nm	317 nm	321 nm
	277 nm	280 nm	279 nm	284 nm
	239 nm	238 nm	237 nm	239 nm
E ^{1%} 1cm	190 376 384	165 346 353	166 382 336	204 260 426
3	16100	14000	14000	17300
	31800	29300	32400	22100
	32500	29900	28400	36100

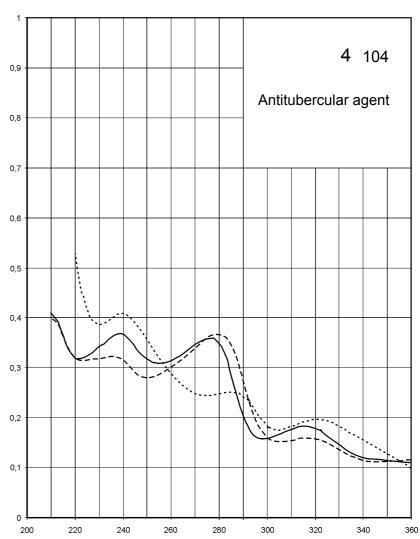


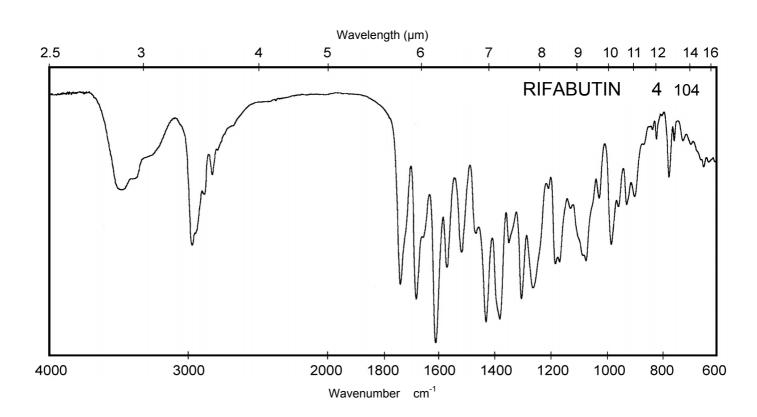

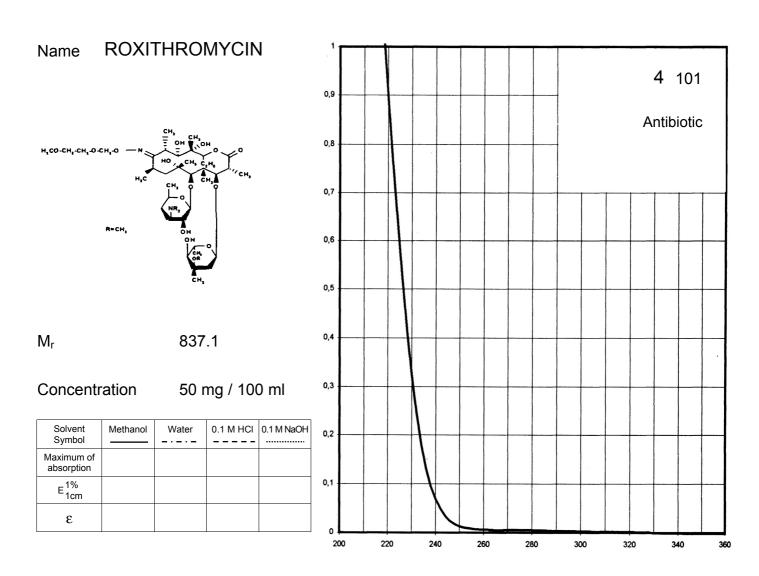

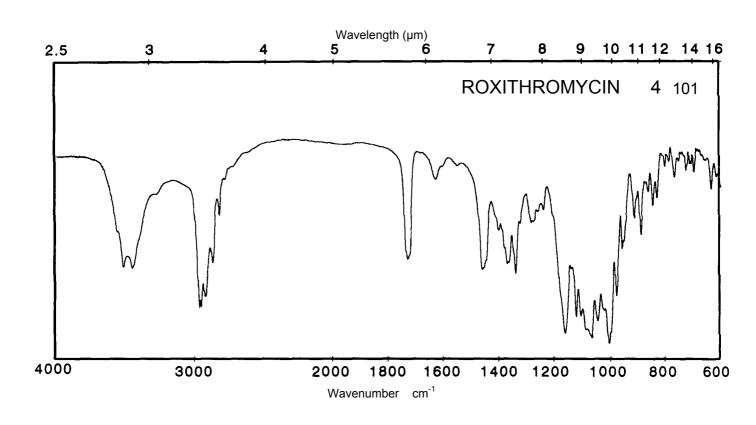



Name RAUBASINE HYDROCHLORIDE

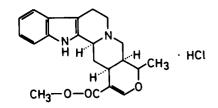

 M_r 388.9

Concentration 0.5 mg / 100 ml 2 mg / 100 ml

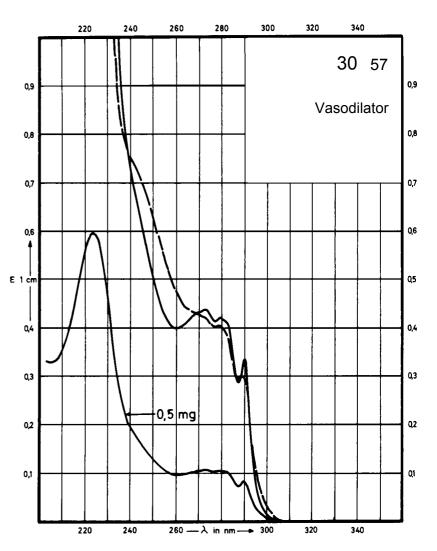

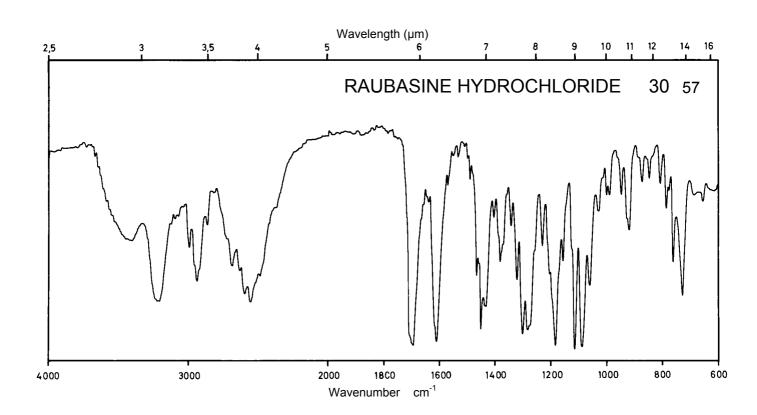



Name RAUBASINE

 M_r 352.4


Concentration 0.5 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm 227 nm		289 nm 278 nm	
E ^{1%} 1cm	227 1260		162 215	
ε	8000 44400		5720 7560	

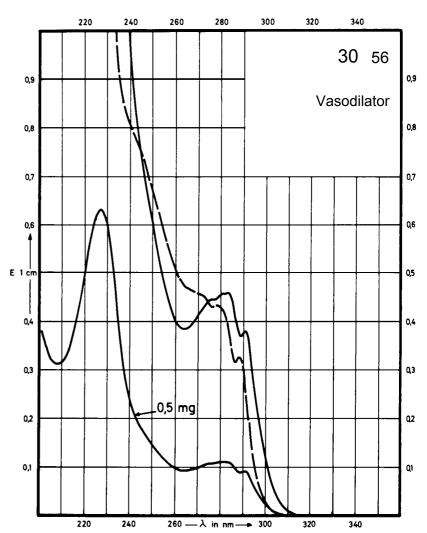

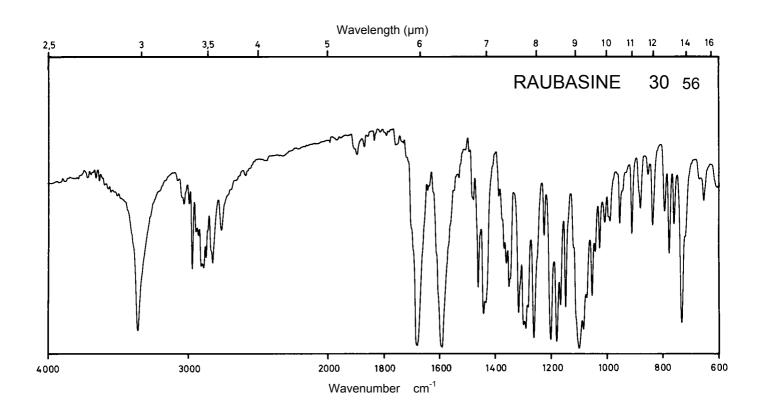



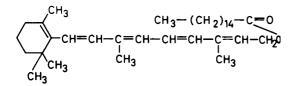



524.9

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	322 nm			
E ^{1%} 1cm	1050			
ε	35590			

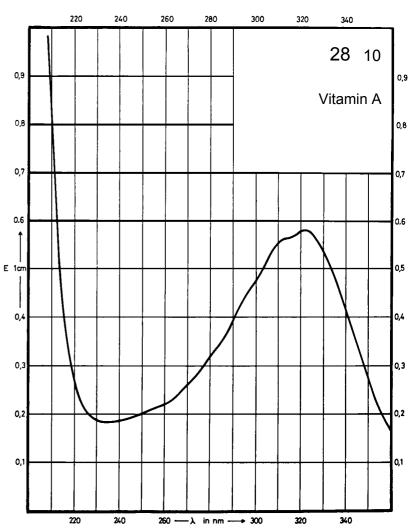

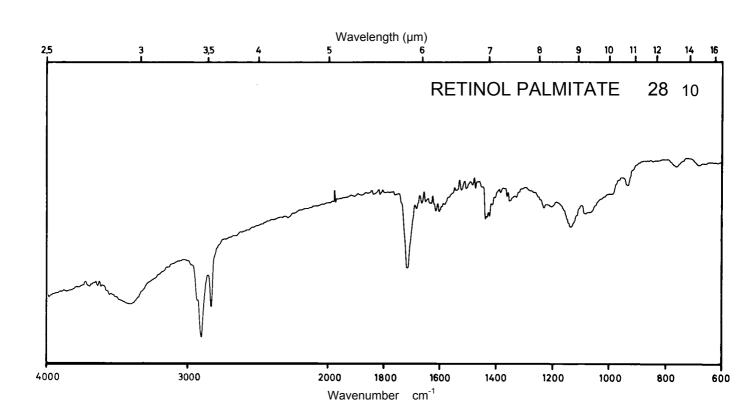




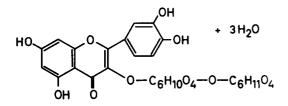


 M_{r}

664.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
-,				
Maximum of absorption	382 nm 264 nm		351 nm 254 nm	330 nm 270 nm
E ^{1%} 1cm	291 378		273 327	219 351
3	19340 25120		18140 21730	14550 23330

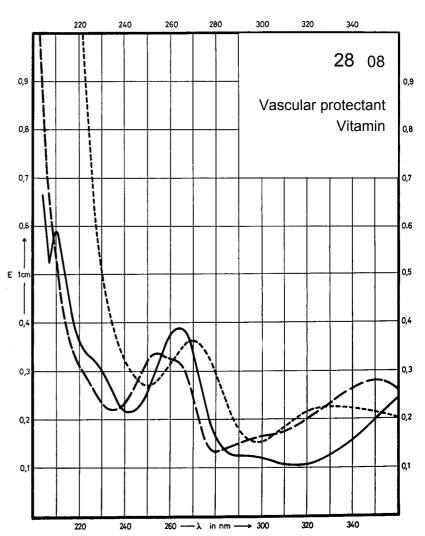

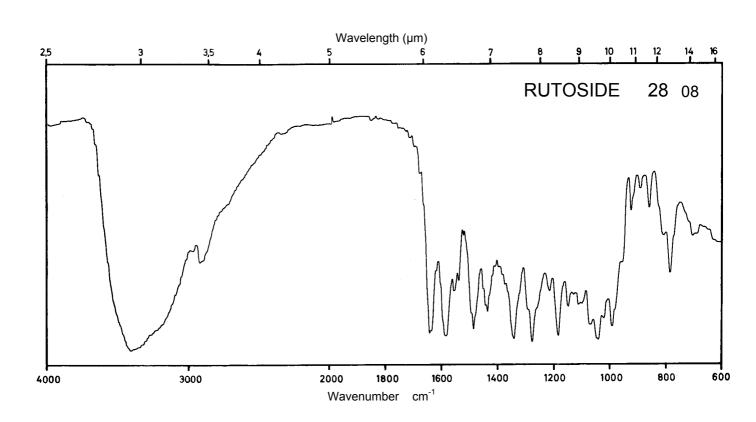

Name RIBOFLAVIN

CH₂OH HO-C-H HO-C-H HO-C-H CH₂ H₃C N N O


 M_{r}


376.4

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	352 nm 269 nm		267 nm 223 nm	356 nm 270 nm
E ^{1%} 1cm	274 890		854 760	325 959
ε	10310 33500		32140 28610	12230 36100

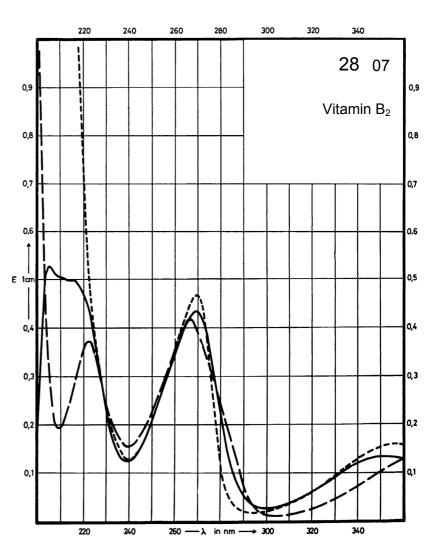

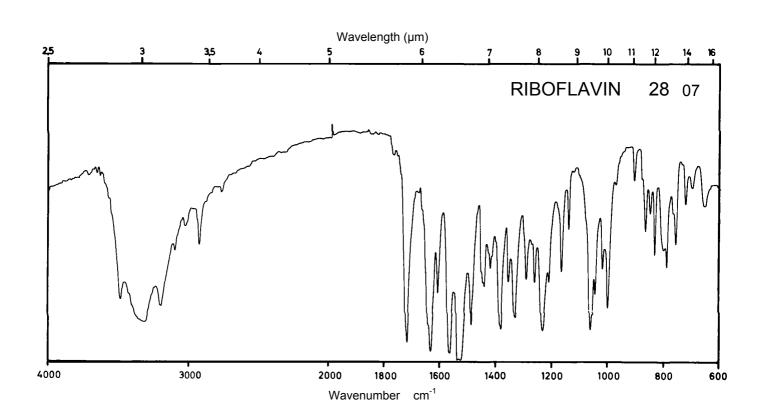





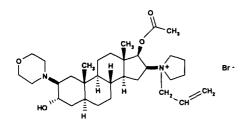


609.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
ε				

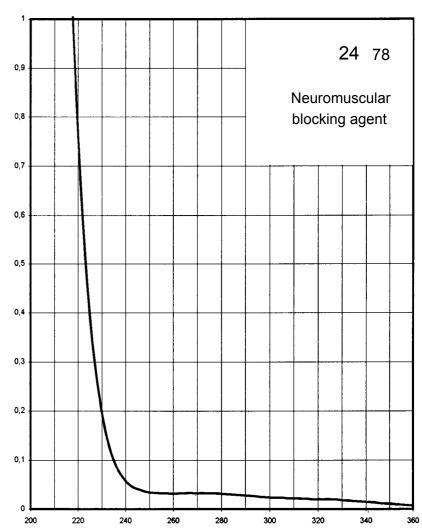

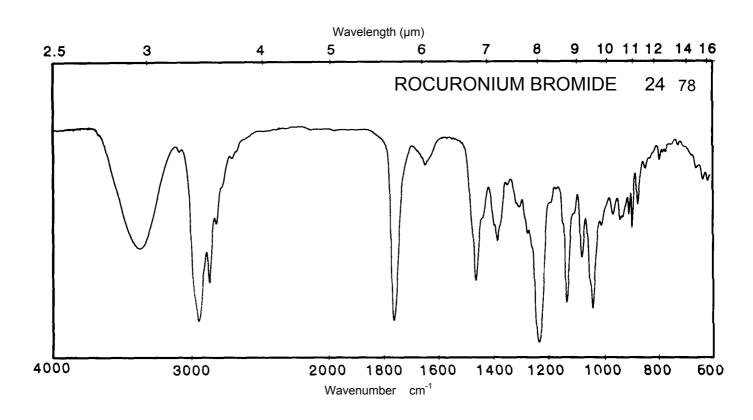



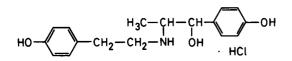



 M_r

323.8


Concentration 1 mg / 100 ml 2 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 277 nm 225 nm 274 nm 222 nm 292 nm 240 nm absorption E^{1%} 1cm 102 565 72 346 3300 18300 2730 15900 2320 11200 ε

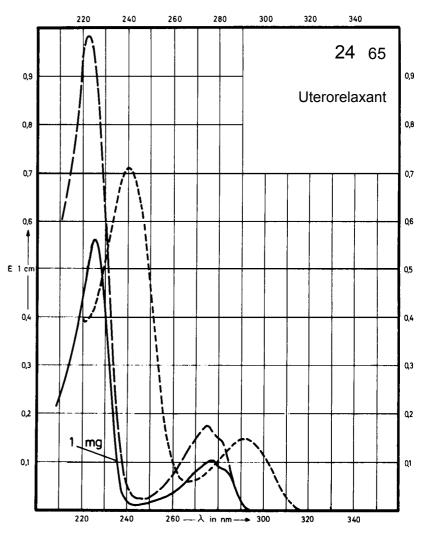

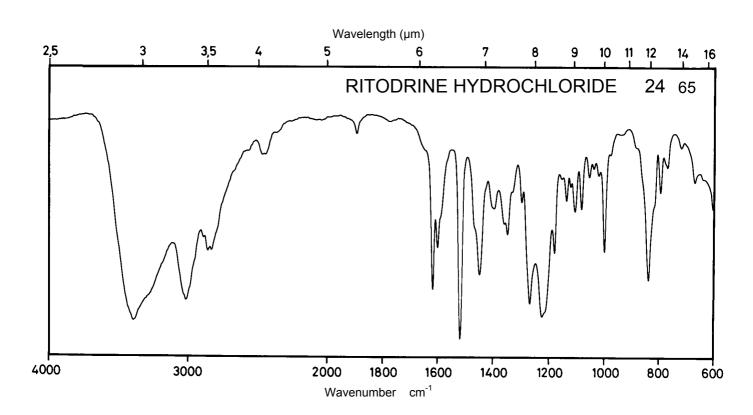





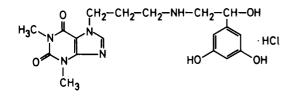


 M_r 425.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm		274 nm	274 nm
E ^{1%} 1cm	250		251	242
ε	10700		10700	10300

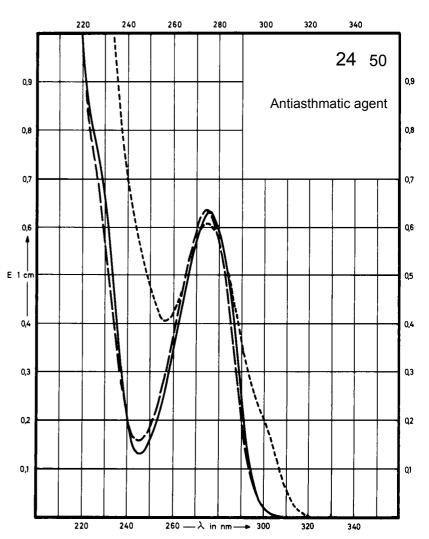

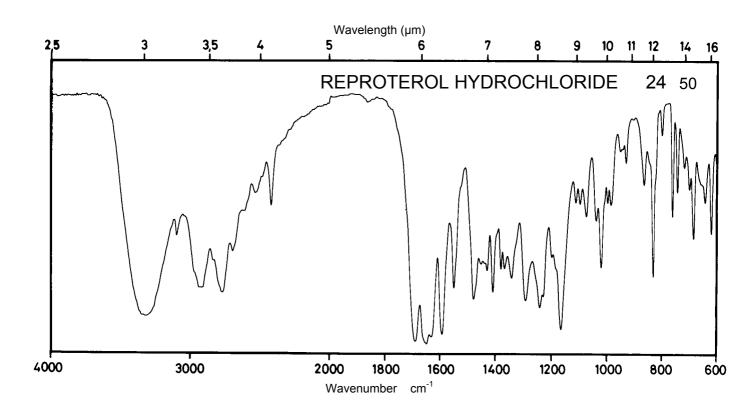




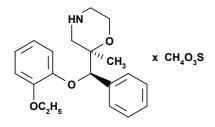


 M_{r} 409.5

Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm	274 nm	274 nm	275 nm
E ^{1%} 1cm	54	52	53	54
ε	2200	2130	2160	2200

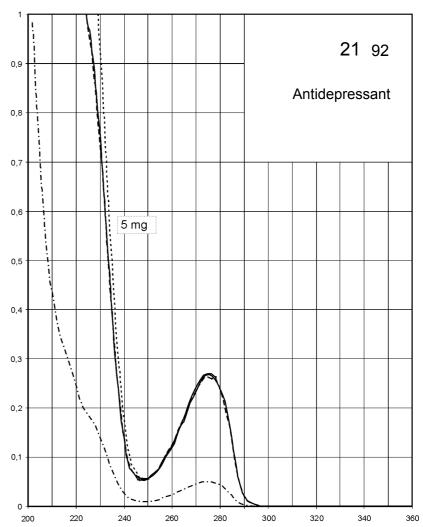

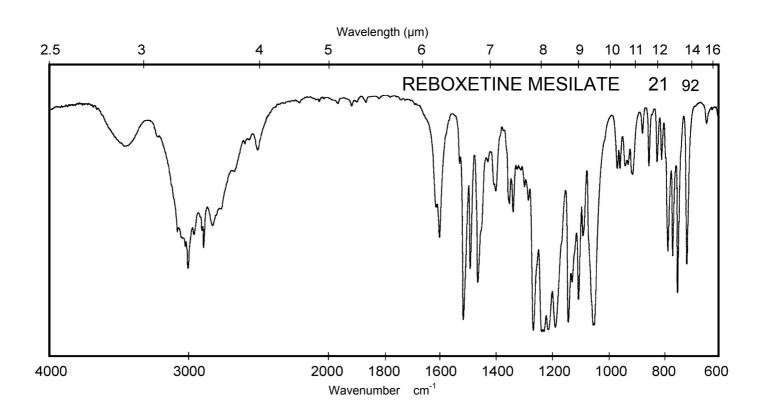



Name ROXATIDINE ACETATE HYDROCHLORIDE

$$\begin{array}{c} \text{OC-CH}_{3} \\ \text{O} \\ \text{CH}_{2} \\ \text{O-(CH}_{2})_{3} - \text{NH-CO} \\ \\ \text{HCI} \end{array}$$


 M_r 384.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 283 nm	275 nm	275 nm	
E ^{1%} 1cm	58 52	50	50	
ε	2220 2000	1930	1930	

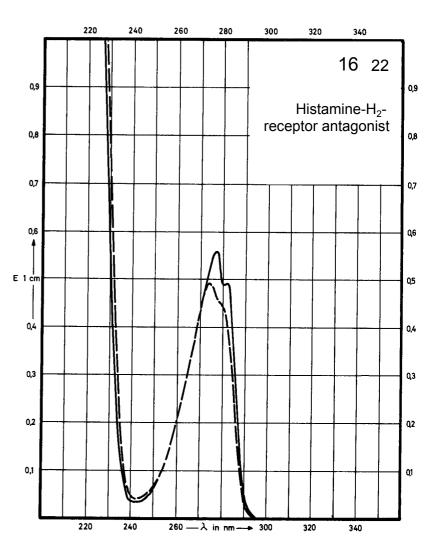

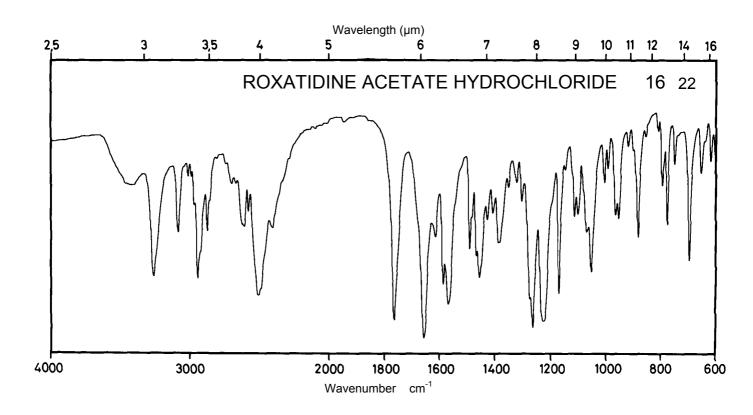



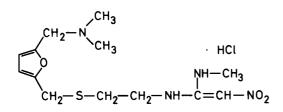



 M_r

Concentration 1 mg / 100 ml

350.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	325 nm 228 nm		317 nm 225 nm	309 nm
E ^{1%} 1cm	527 490		31 587	369
ε	18500 17200		1100 20600	12900

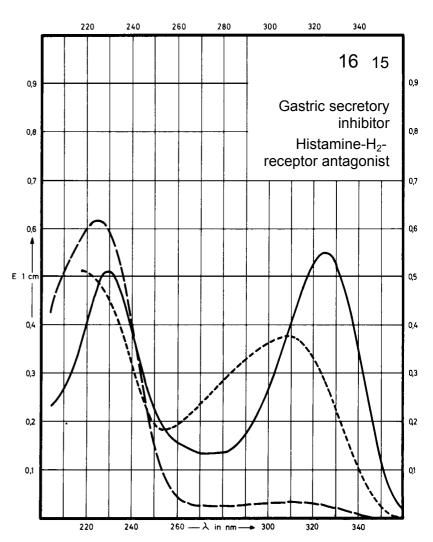

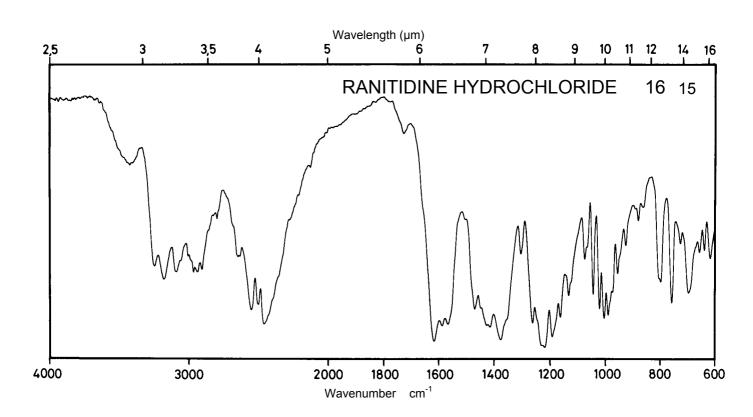








 M_{r}


110.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm		272 nm	289 nm
E ^{1%} 1cm	194		169	312
ε	2140		1860	3440

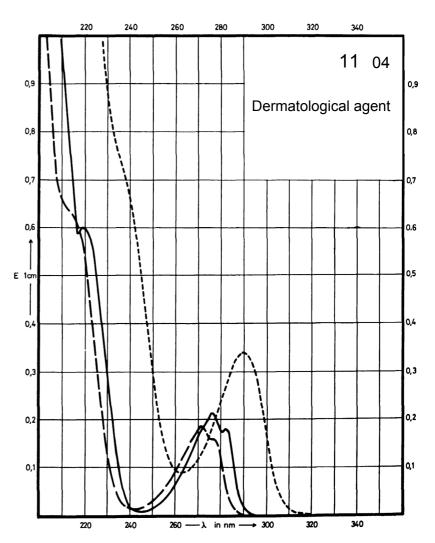


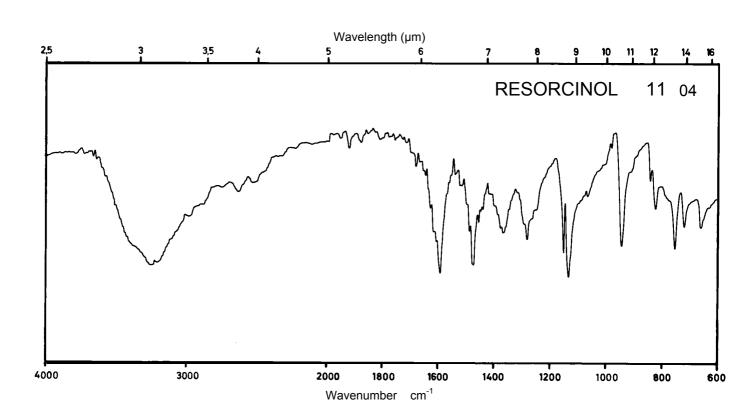

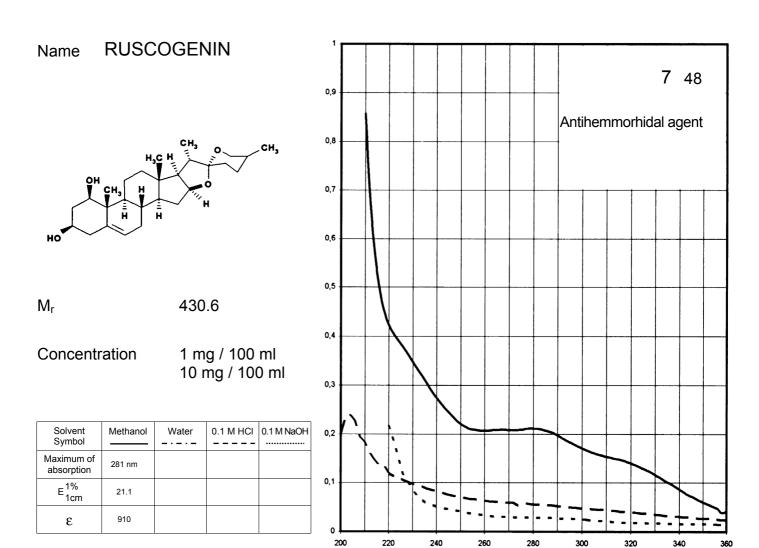

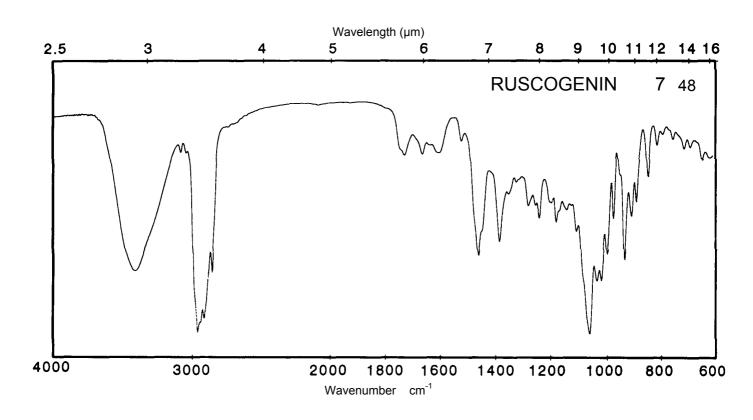




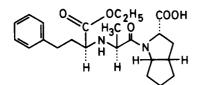


416.5


80 mg / 100 ml Concentration

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm		258 nm	258 nm
E ^{1%} 1cm	5.0		4.6	5.1
ε	210		190	213

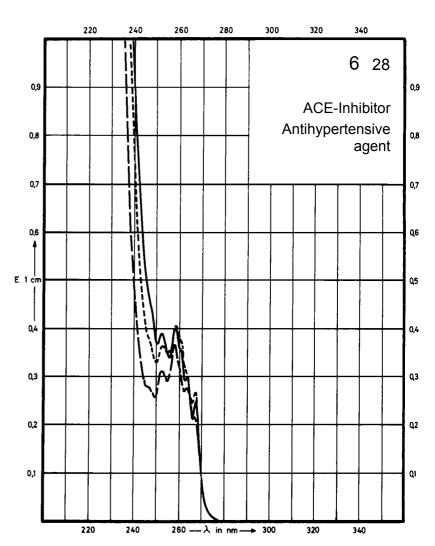

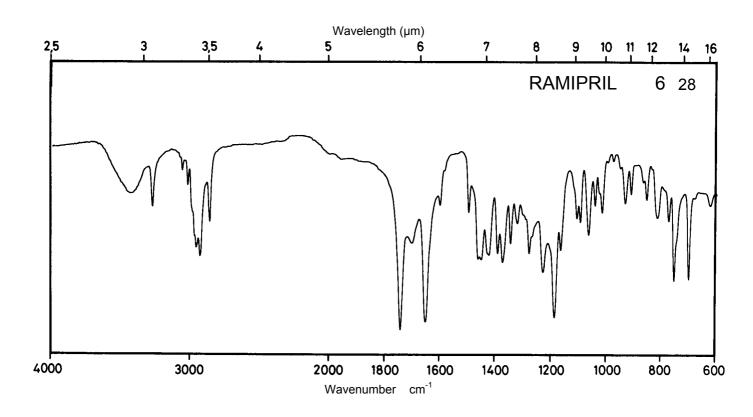



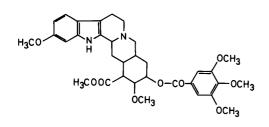


 M_r


608.7

Concentration

2.4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	294 nm 266 nm		290 nm 268 nm	
E ^{1%} 1cm	168 274		178 278	
ε	10230 16680		10830 16920	

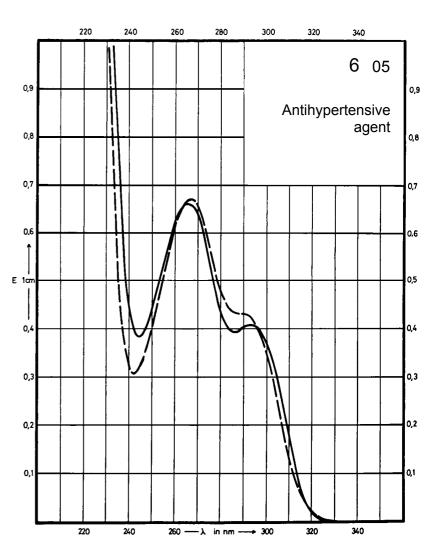

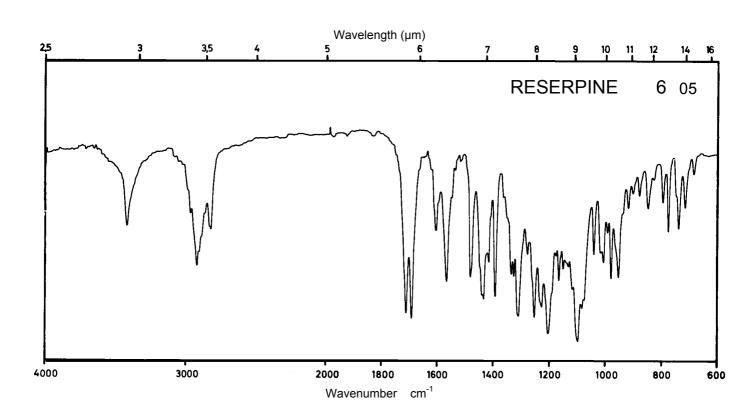

Name ROSOXACIN


 M_{r}


294.3

Concentration

0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	317 nm 269 nm		277 nm	337 nm 270 nm
E ^{1%} 1cm	329 1600		1580	294 1660
ε	9700 47000		46500	8600 48800

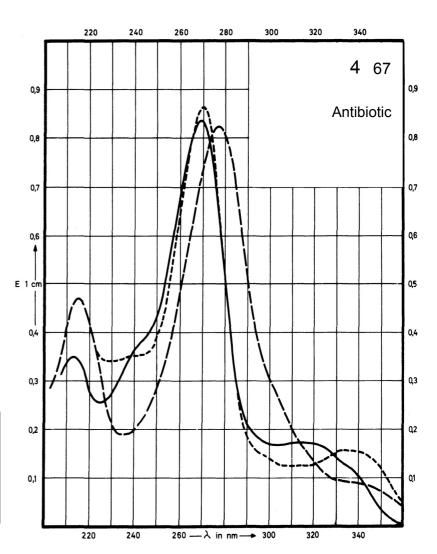

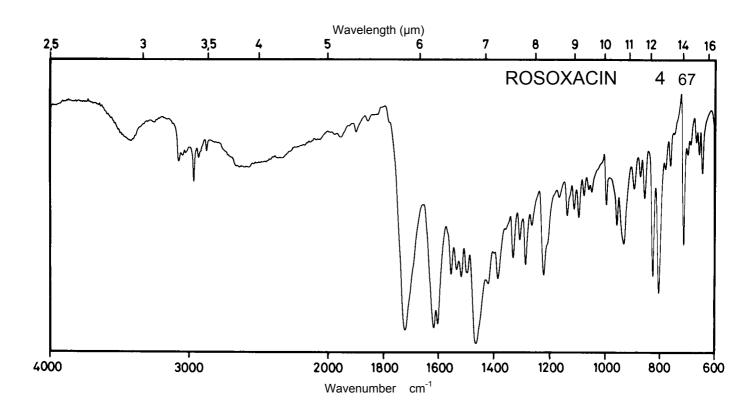



Name RIFAMPICIN

 M_{r}


823.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	480 nm 337 nm 237 nm		470 nm 336 nm 231 nm	470 nm 332 nm
E ^{1%} 1cm	190 340 420		116 250 320	176 310
3	15640 27980 34570		9550 20520 26340	14480 25510

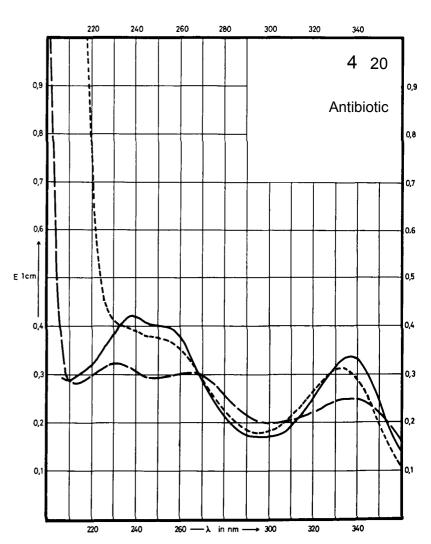

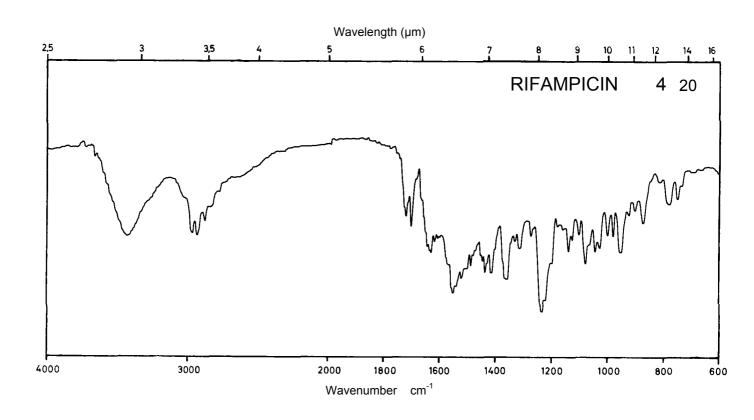

Name ROLITETRACYCLINE

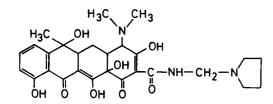


 M_{r}


527.6

Concentration


2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm 220 nm		356 nm 272 nm	268 nm 248 nm
E ^{1%} 1cm	349 265		288 339	266 261
3	18410 13980		15190 17890	14030 13770

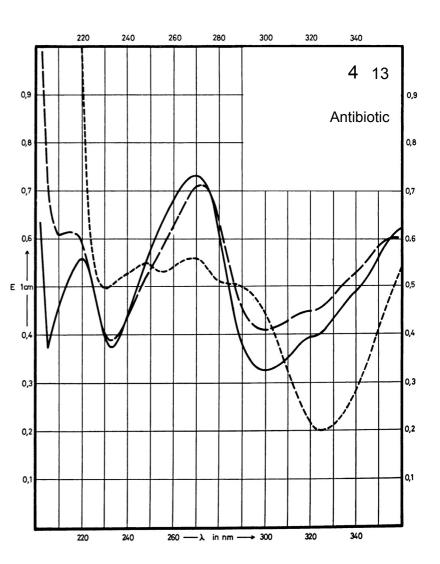

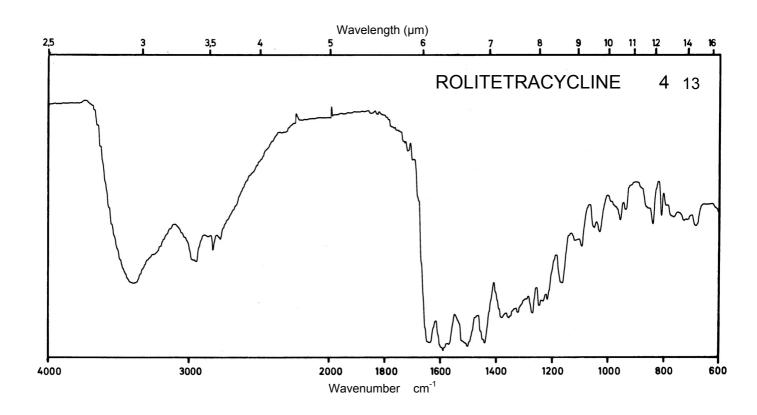




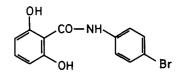


 M_{r}


308.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm			289 nm
E ^{1%} 1cm	750			686
ε	23100			21100

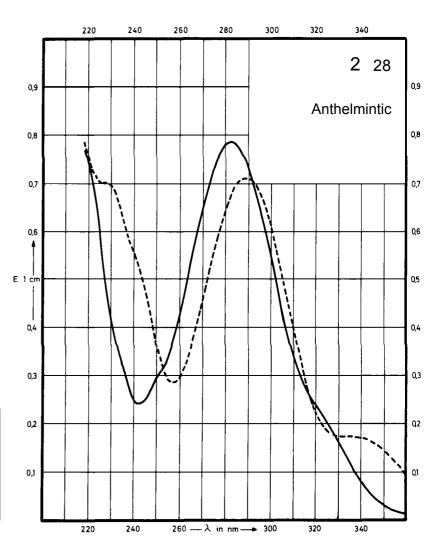

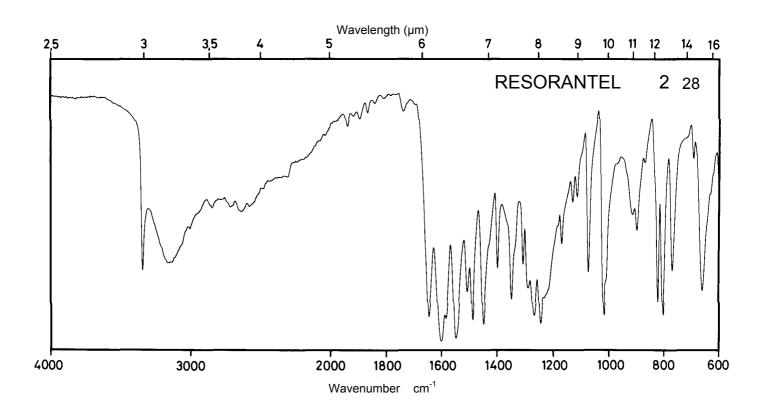




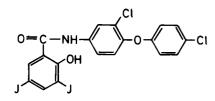


 M_{r}


626.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	363 nm 281 nm 239 nm			362 nm 284 nm 235 nm
E ^{1%} 1cm	139 306 524			158 283 564
3	8700 19160 32800			9890 17720 35310

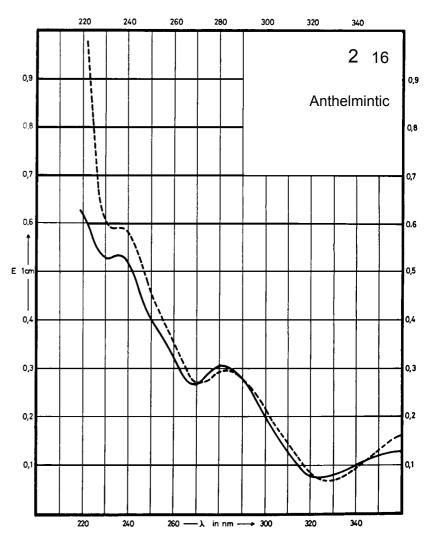

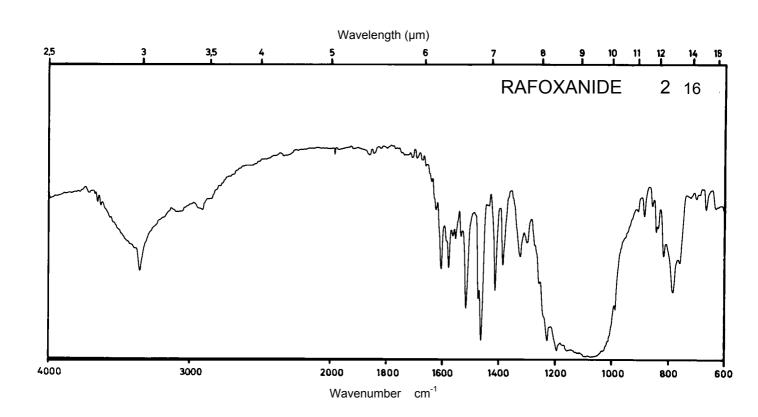




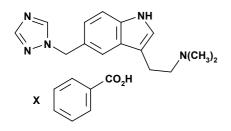


391.5


Concentration 0.5 mg/100 ml 2.5 mg/100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm	278 nm	281 nm	281 nm
	227 nm	225 nm	226 nm	226 nm
E ^{1%}	148	139	152	129
1cm	1293	1227	1240	1204
ε	5780	5430	5930	5060
	50600	48100	48600	47100

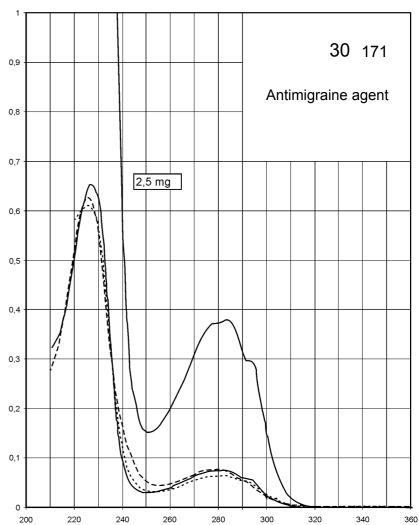

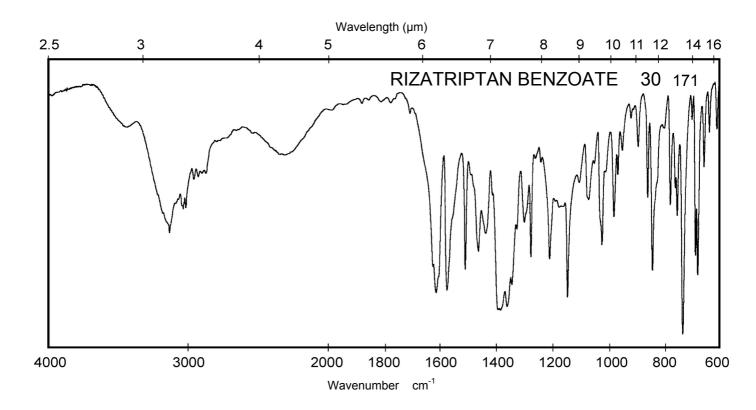





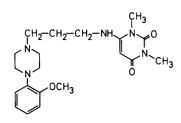



387.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm	269 nm 234 nm	268 nm	269 nm 235 nm
E ^{1%} 1cm	664	669 237	652	683 243
ε	25730	25930 9200	25280	26440 9410

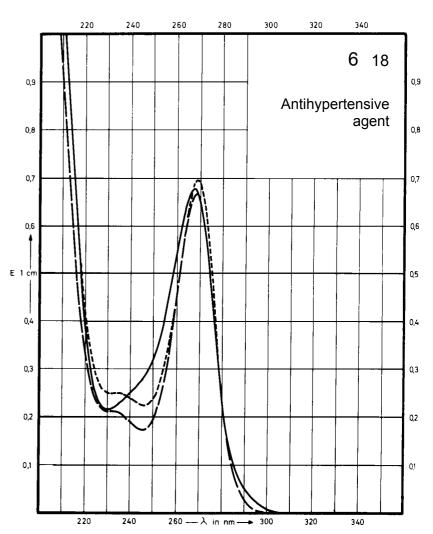

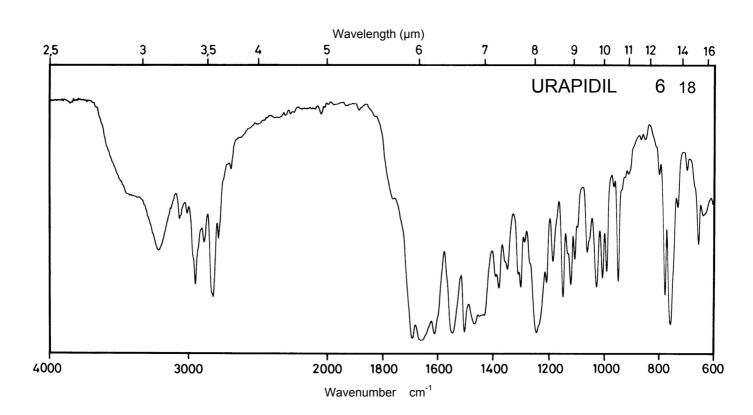

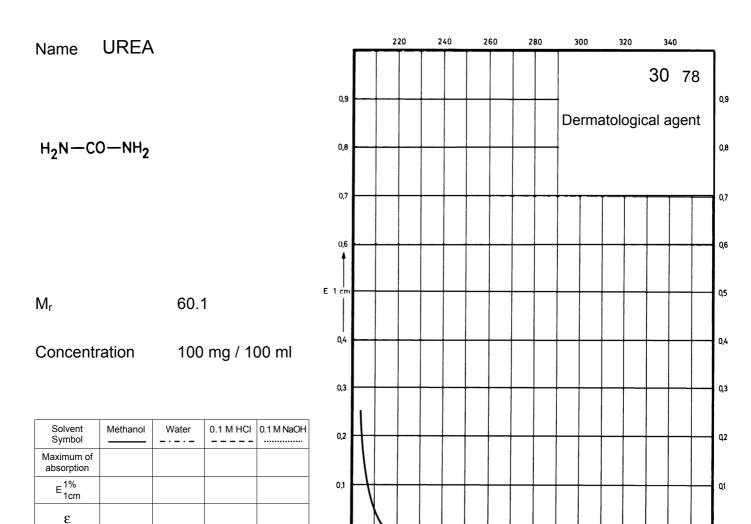

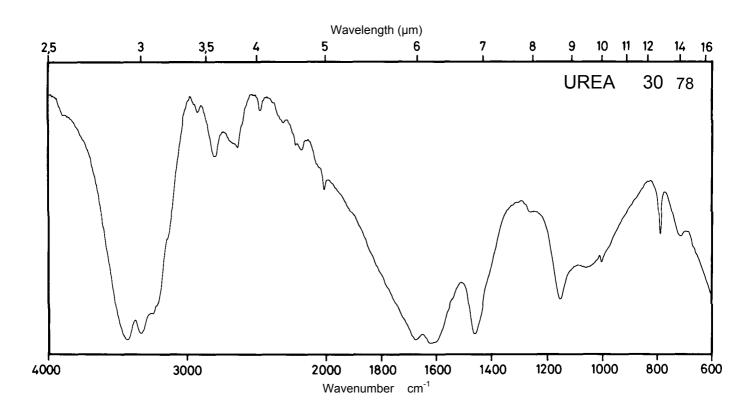

260 — λ in nm → 300

340


Name URIC ACID


HO N N OF

 M_{r}


168.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	293 nm 235 nm		283 nm 228 nm	295 nm
E ^{1%} 1cm	740 570		720 530	800
3	12440 9580		12100 8910	13450

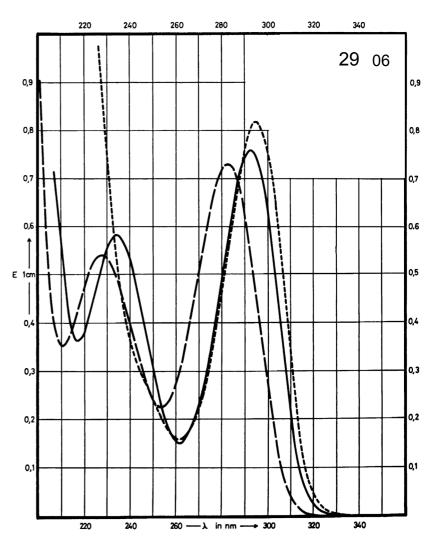

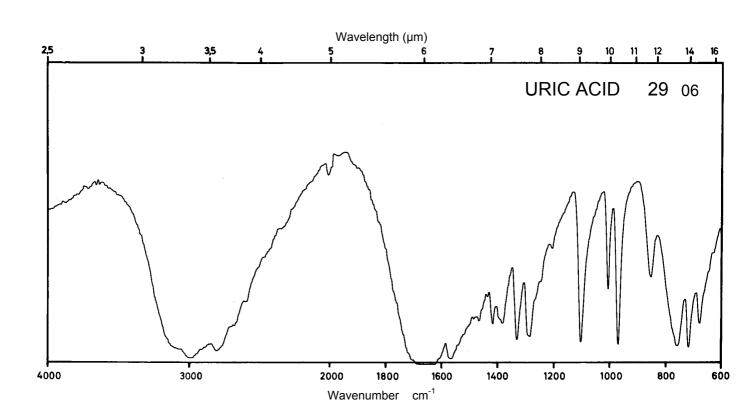




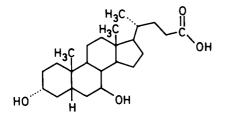


 M_r

392.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

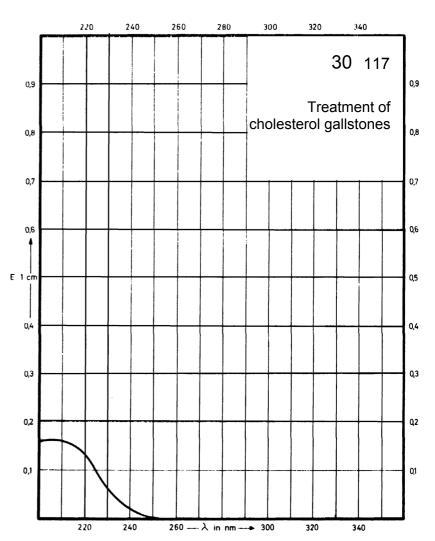

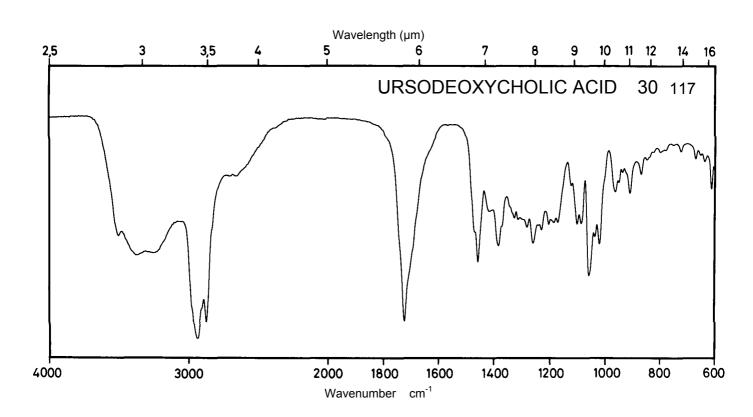


Name VINCAMINE


CH3O-CO OH CH2-CH3

 M_{r}


354.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm		268 nm	282 nm
E ^{1%} 1cm	239		242	208
ε	8470		8580	7370

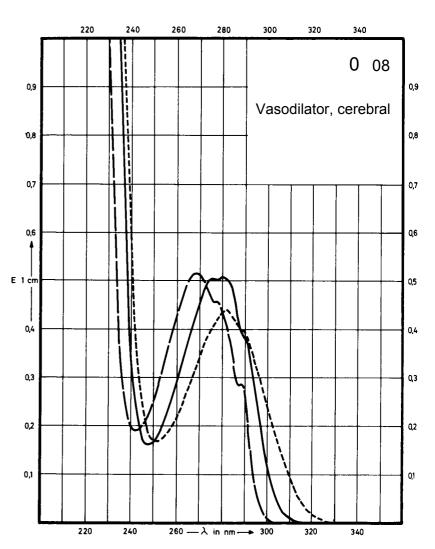

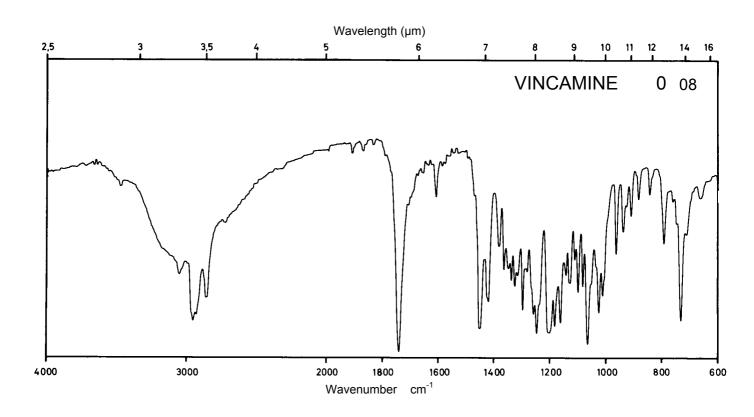





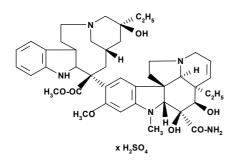


 M_r 852.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 216 nm	271 nm 215 nm	270 nm 215 nm	268 nm
E ^{1%} 1cm	38.3 115	35.6 113	36.6 118	31.4
ε	3260 9830	3030 9660	3030 10000	2680

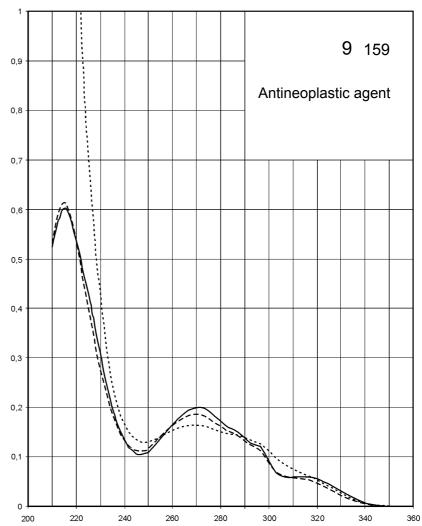

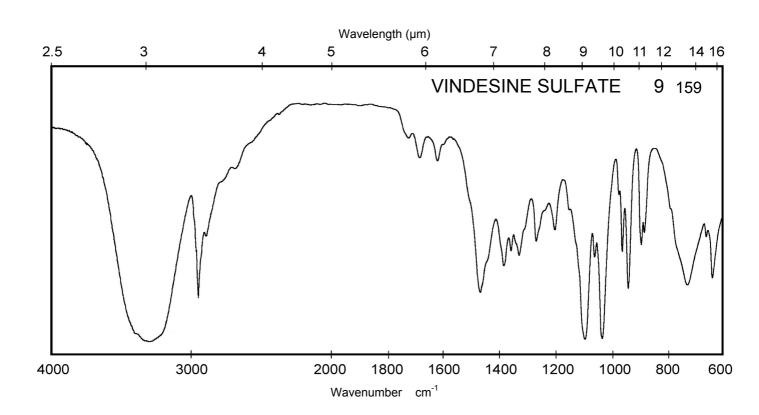





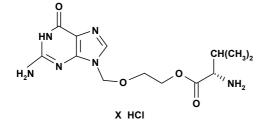


 M_r 360.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	255 nm	252 nm	256 nm	259 nm
E ^{1%} 1cm	404	375	346	310
ε	14600	13500	12500	11200

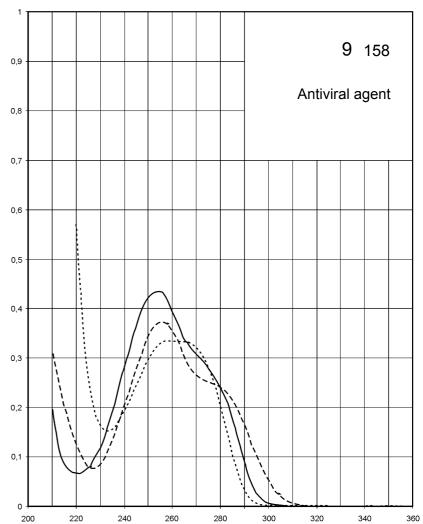

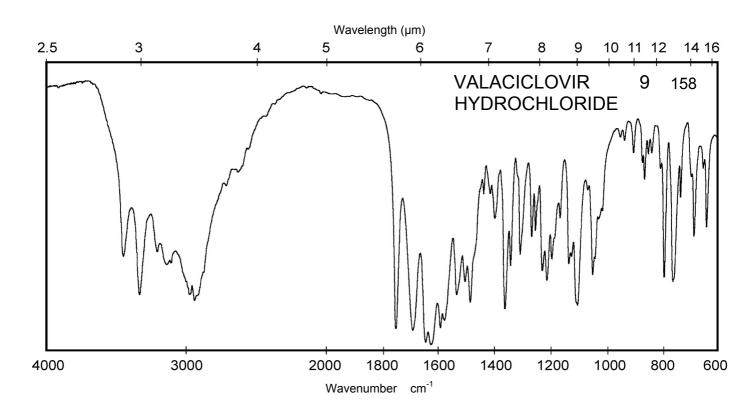

Name VIDARABINE

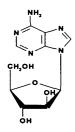


 M_{r}


267.2

Concentration

1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm	259 nm	257 nm	259 nm
E ^{1%} _{1cm}	574	570	558	579
ε	15300	15200	14900	15500

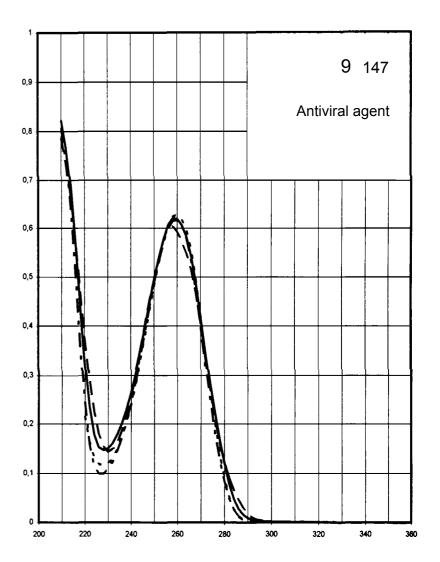

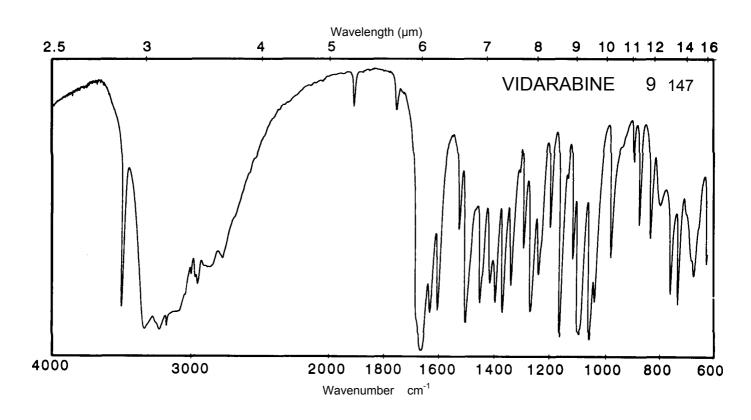







152.1


Concentration 0.2 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 307 nm		278 nm 308 nm	248 nm 348 nm
E ^{1%} 1cm	685 660		685 610	610 1640
ε	10420 10040		10420 9280	9280 24940

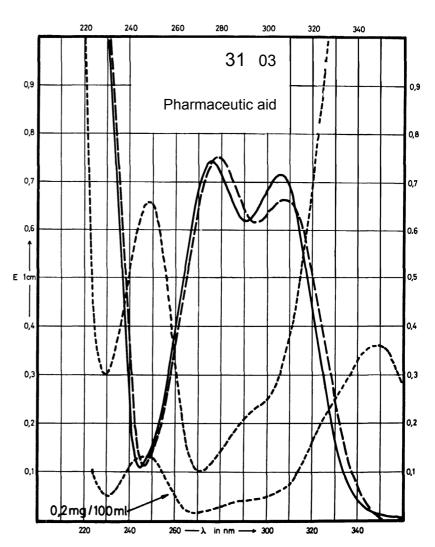

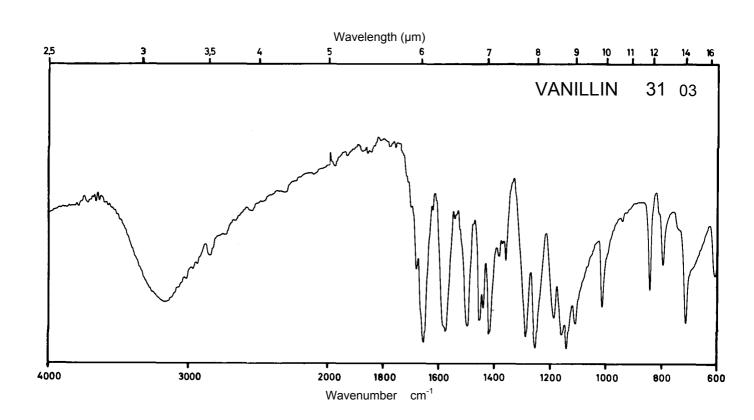





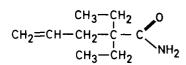



155.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

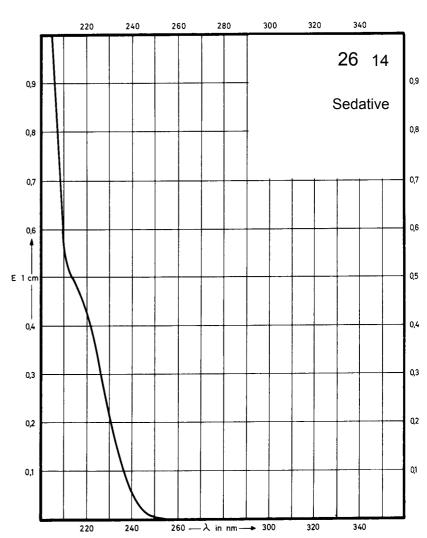

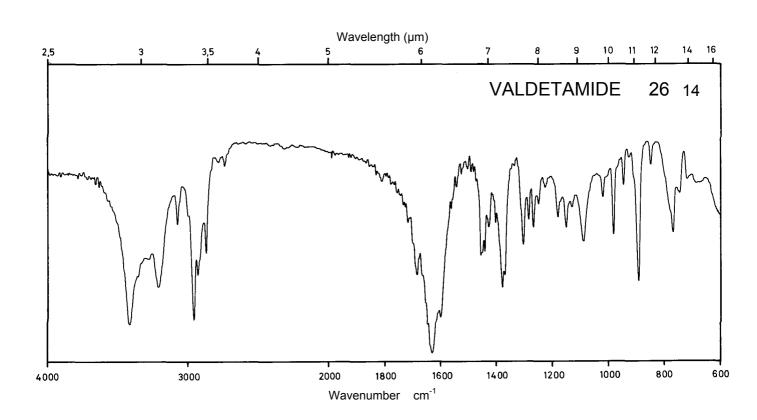


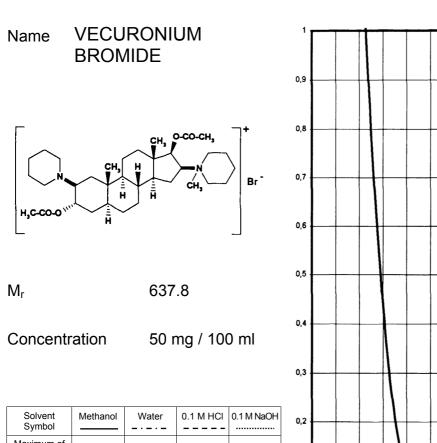

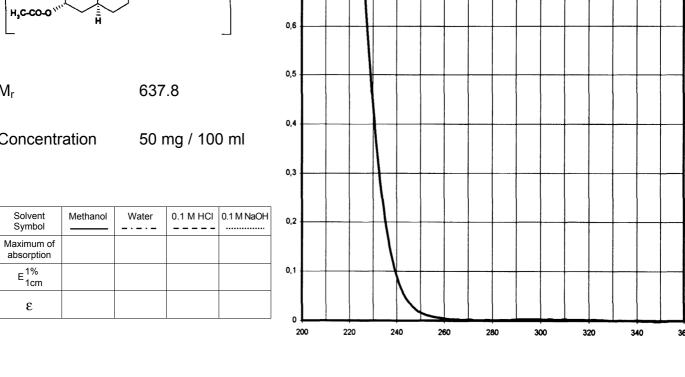


24 81

Neuromuscular blocking agent

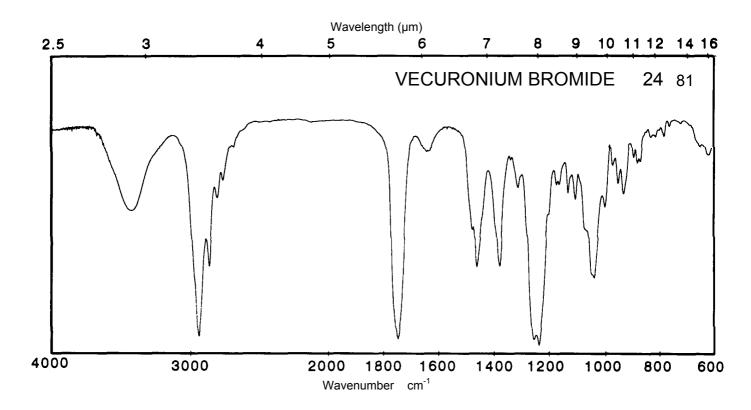




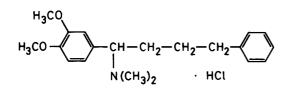


 M_{r}


349.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	279 nm 235 nm		277 nm 233 nm	277 nm
E ^{1%} 1cm	90 261		85 259	80
ε	3160 9130		2970 9050	2800

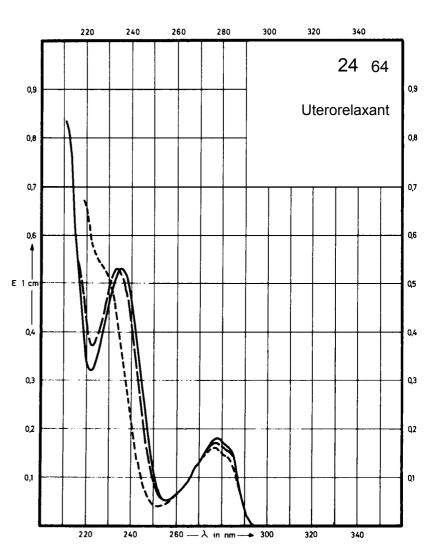

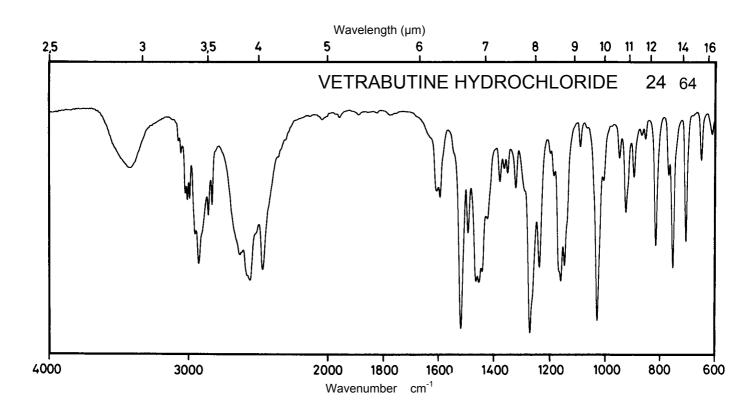




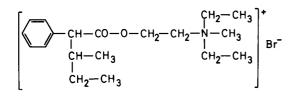


 M_{r} 386.4


Concentration 10 mg / 100 ml 100 mg / 100 ml

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	265 nm		265 nm	
	259 nm		259 nm	
absorption	253 nm		253 nm	
1%	4.10		4.20	
E ^{1%} 1cm	5.26		5.50	
1cm	4.21		4.50	
	160		165	
3	200		215	
_	165		175	

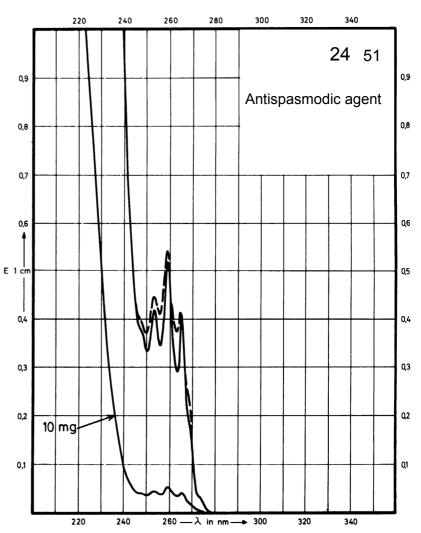

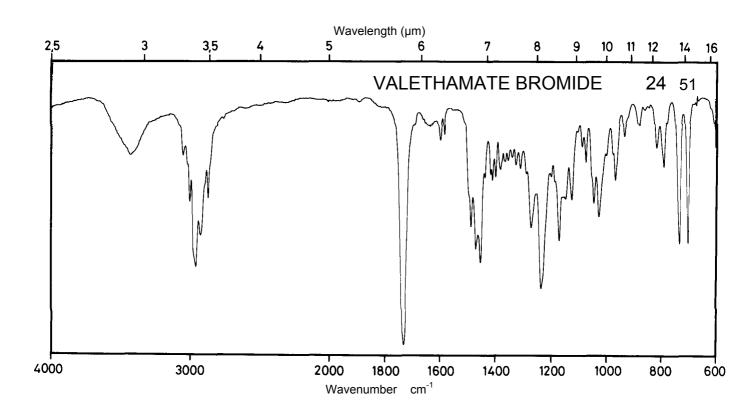


Name VILOXAZINE HYDROCHLORIDE

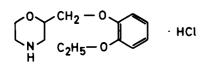



 M_r 273.8


Concentration 2 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm 222 nm		273 nm 221 nm	273 nm
E ^{1%} 1cm	81 270		80 267	82
3	2220 7400		2190 7300	2250

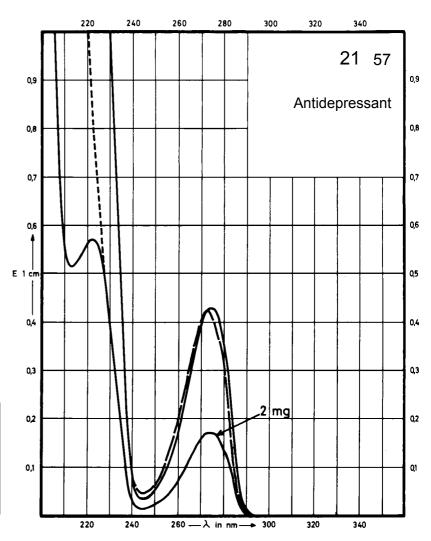


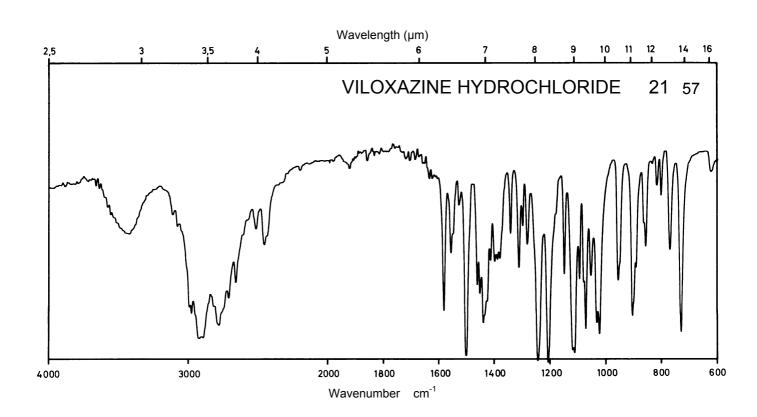

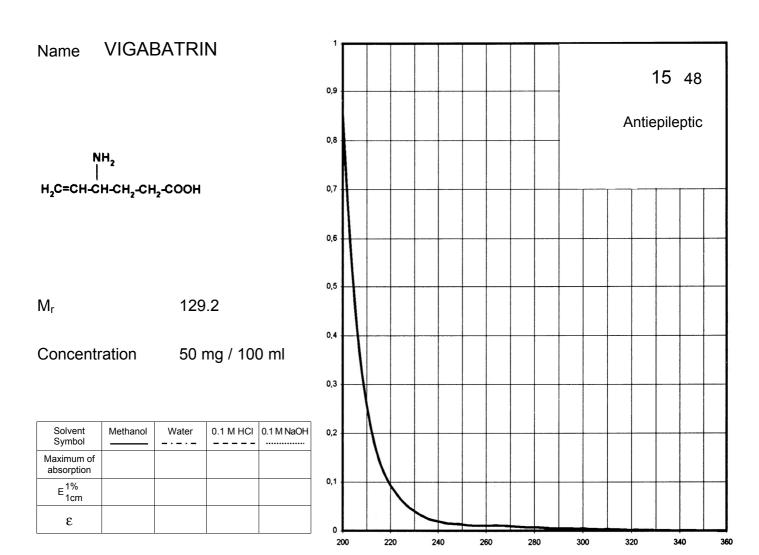

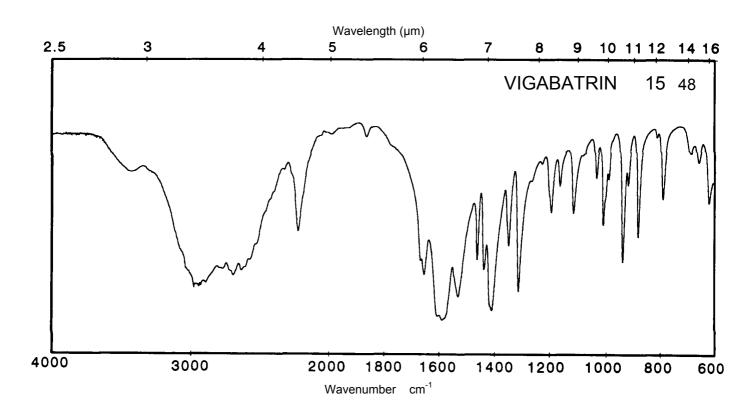





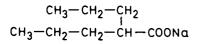

 M_{r} 166.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

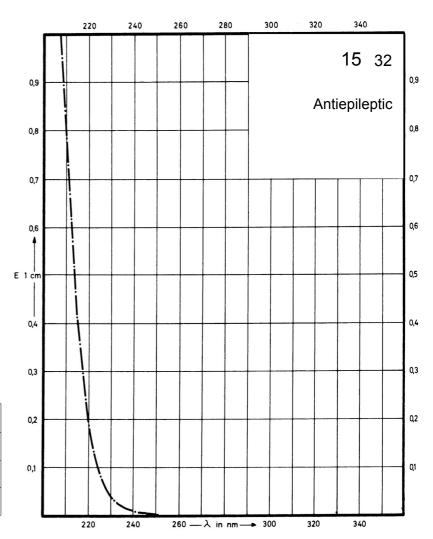

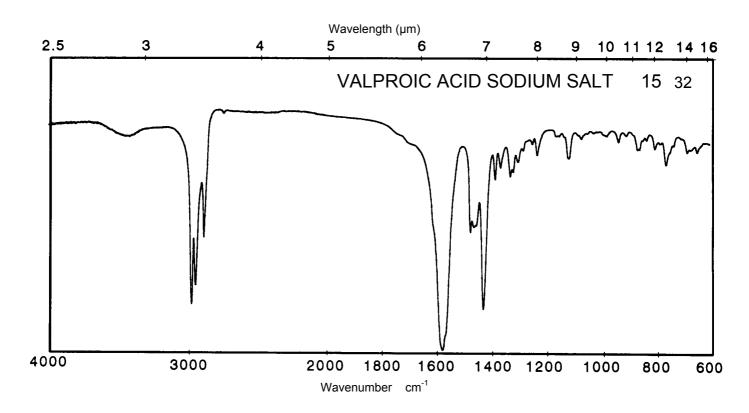


Name VINYLBITAL


$$CH_2 = CH$$
 , $CO - NH$, $CH_3 - CH_2 - CH_2 - CH$, $CO - NH$, CH_3

 M_{r}


224.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				247 nm
E ^{1%} 1cm				298
3				6680

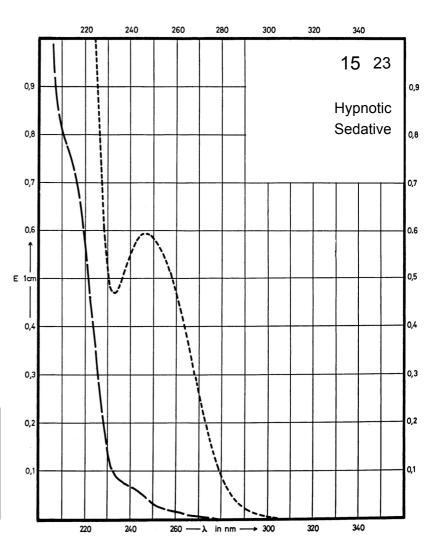

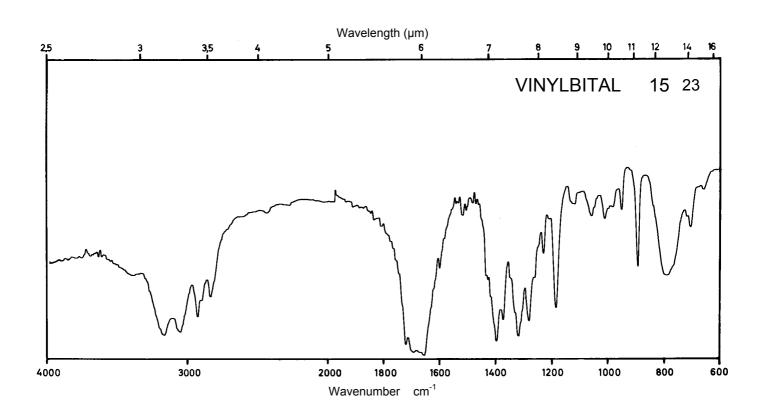





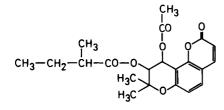


 M_r 388.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	324 nm		324 nm	Decom- position observed
E ^{1%} 1cm	357		355	
ε	13870		13780	

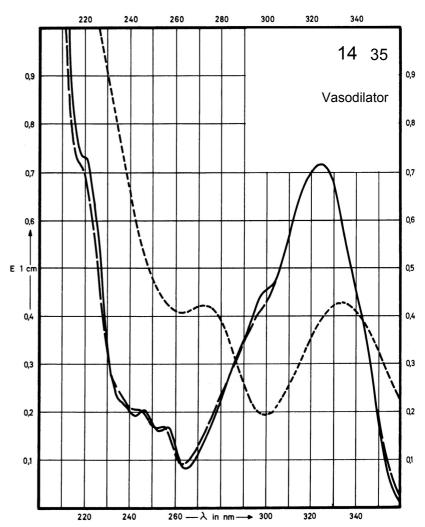

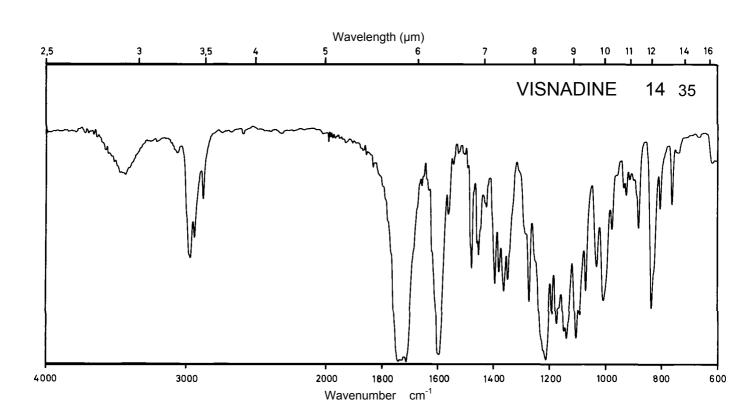


Name VERAPAMIL HYDROCHLORIDE

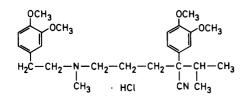



 M_{r}


491.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	279 nm 230 nm		278 nm 228 nm	278 nm
E ^{1%} 1cm	123 335		118 313	113
ε	6050 16400		5800 15400	5550

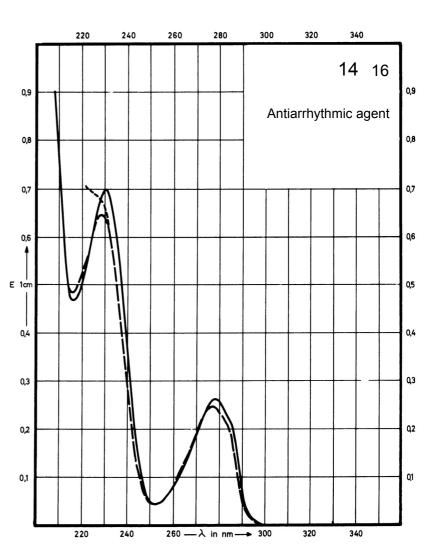


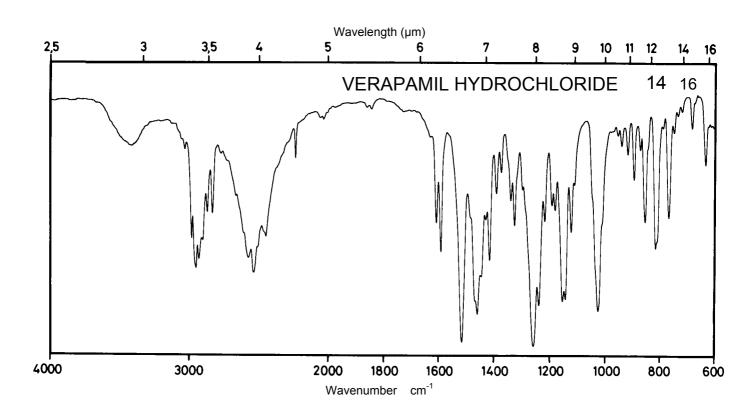

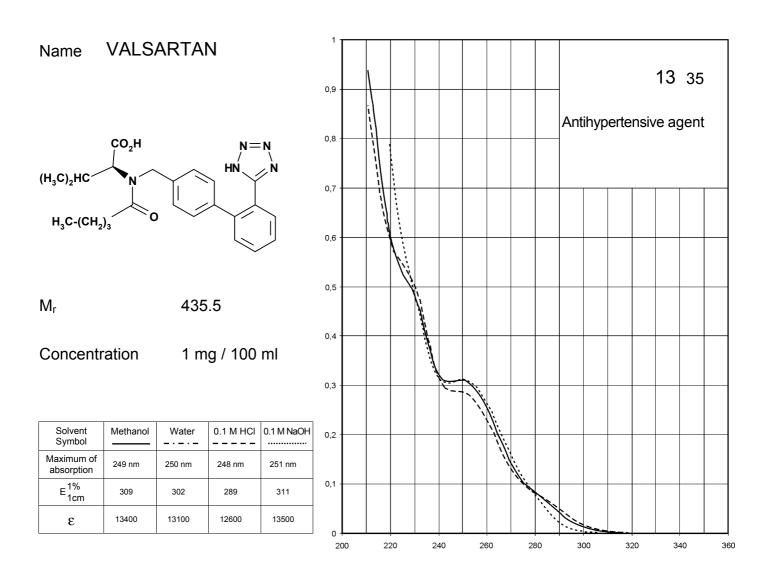

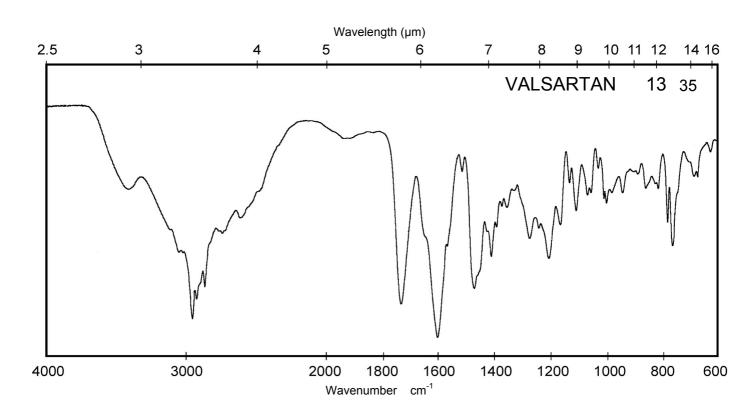

Name VINCRISTINE SULFATE

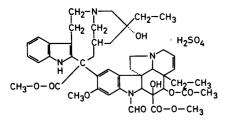


 M_{r}


923.1

Concentration


1.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	297 nm 257 nm 221 nm		297 nm 256 nm 221 nm	299 nm 256 nm
E ^{1%} 1cm	160 161 504		160 162 520	165 163
ε	14810 14870 46490		14810 14970 47980	15250 15050

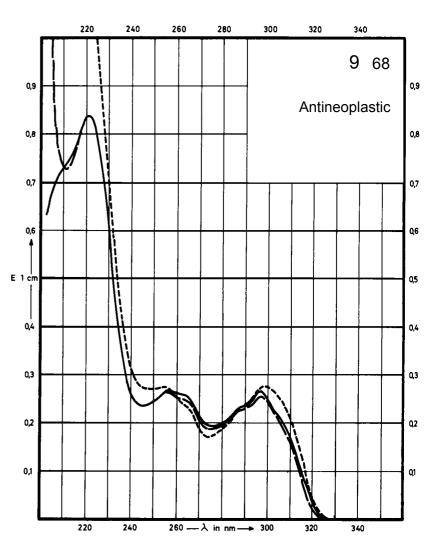

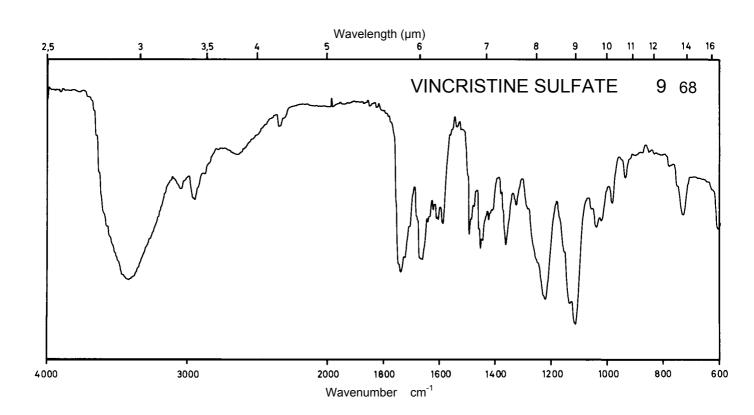



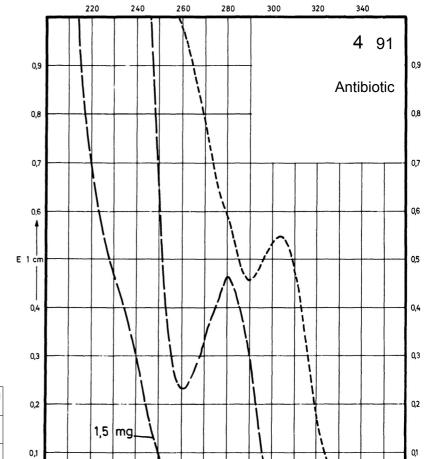


260 — λ in nm →

 M_r 1449.2

Concentration 1.5 mg / 100 ml 10 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		280 nm	280 nm	304 nm
E ^{1%} 1cm		46	47	54
ε		6660	6860	7830

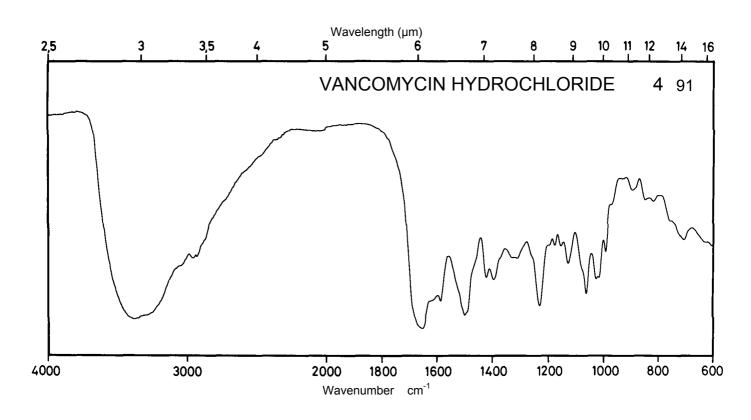

220

240




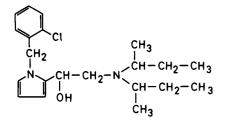


 M_r 362.9


Concentration 2 mg / 100 ml 13 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
ε				

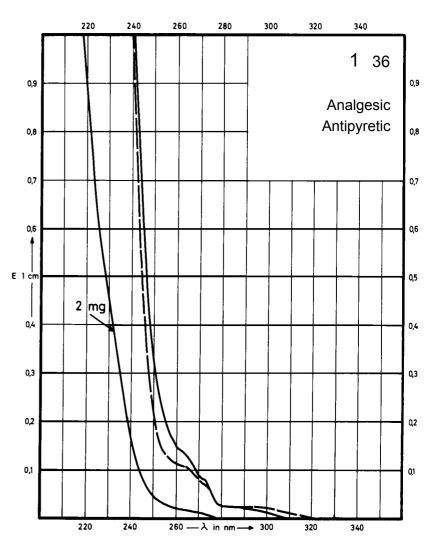

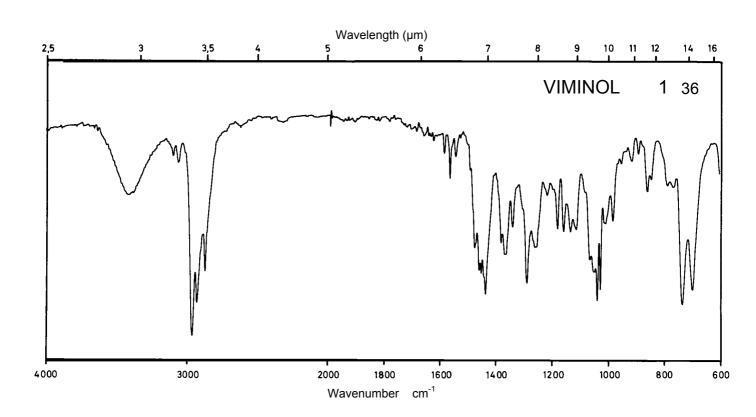



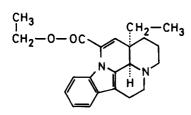

 M_r

Concentration


ε


11800 31500

350.5


1 mg / 100 ml

11350 31300

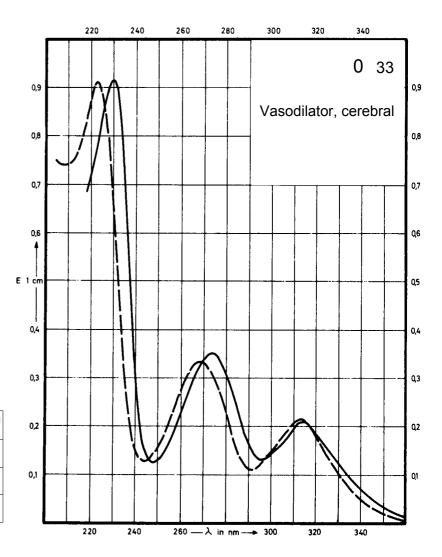

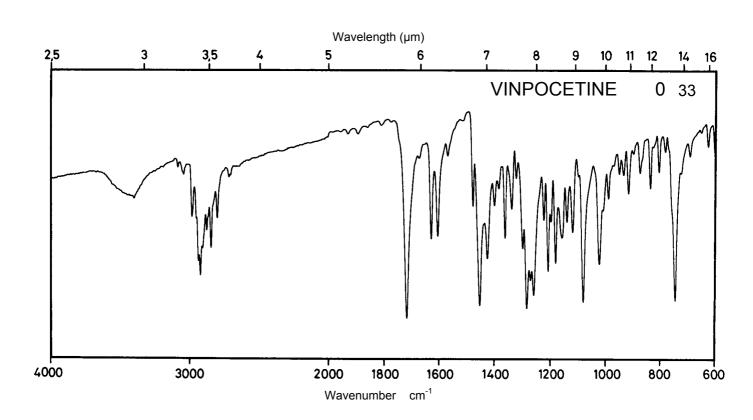


Name VIQUIDIL HYDROCHLORIDE

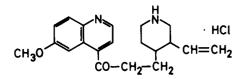


 M_{r} 360.9


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	346 nm	342 nm	355 nm 255 nm	340 nm
E ^{1%} 1cm	129	127	122 508	125
ε	4650	4580	4400 18300	4530

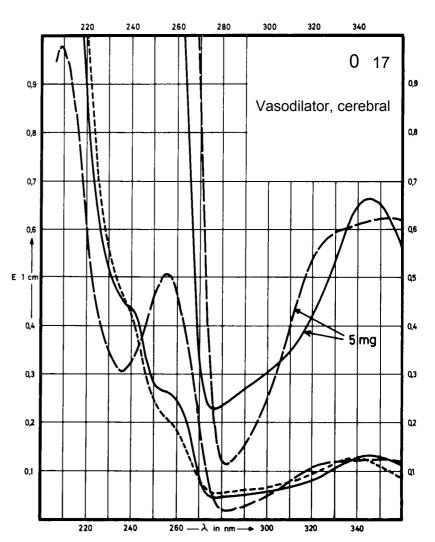

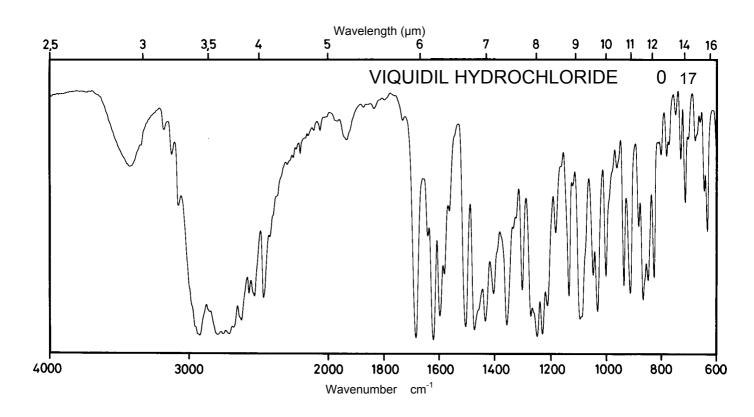




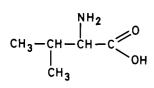


 $\,M_r\,$

117.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

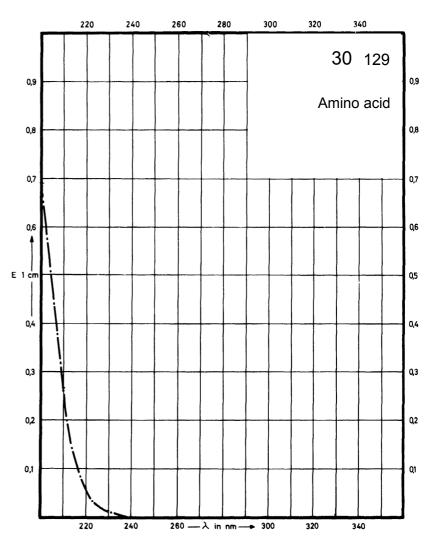

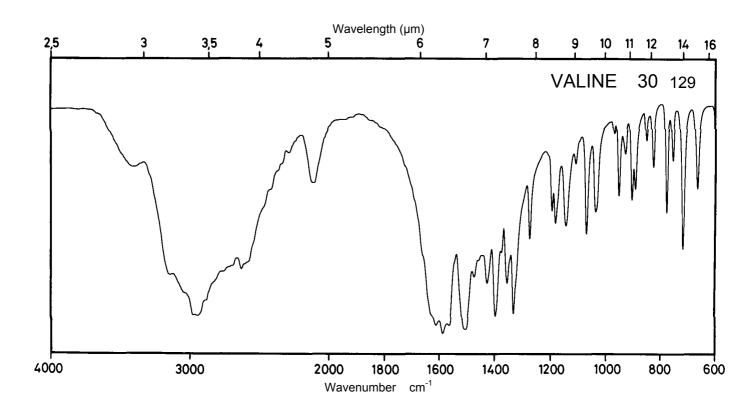

Name XIPAMIDE

H₂NO₂S OH CH₃


 M_{r}


354.8


Concentration

0.8 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	318 nm		296 nm	330 nm
E ^{1%} 1cm	100		92	160
ε	3550		3260	5680

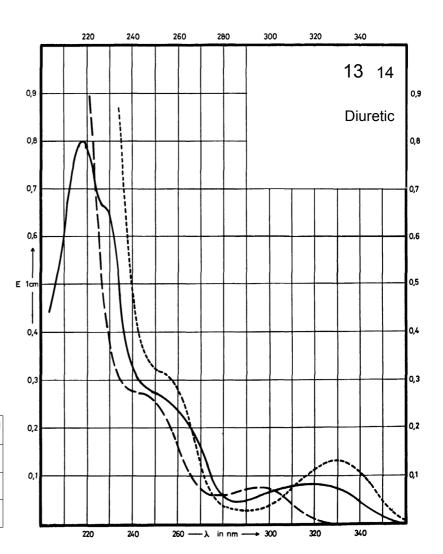


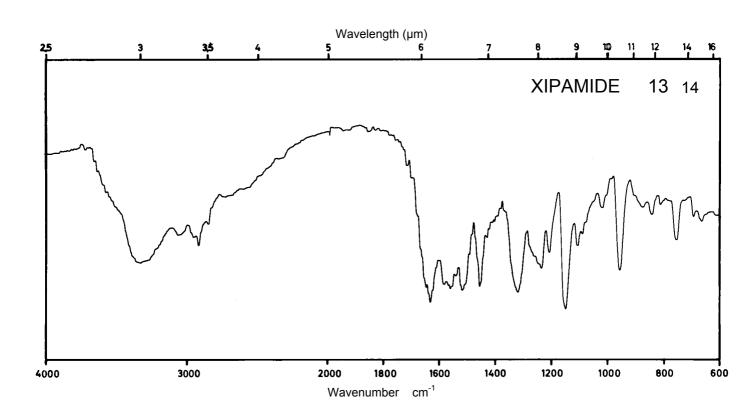

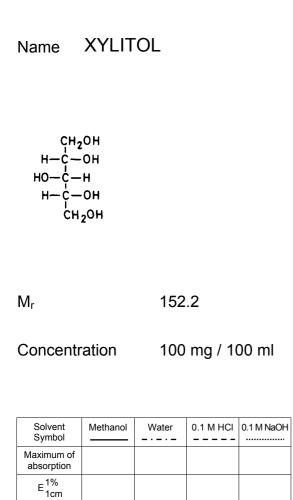


3

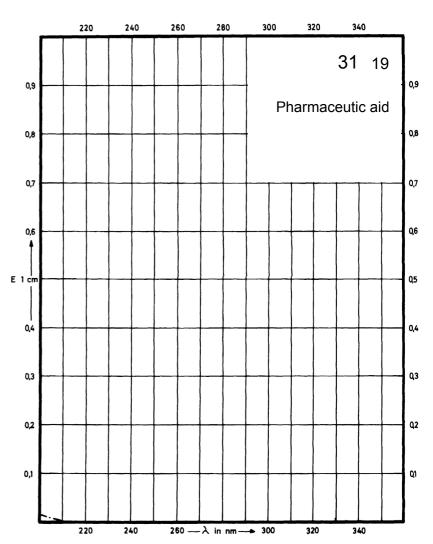

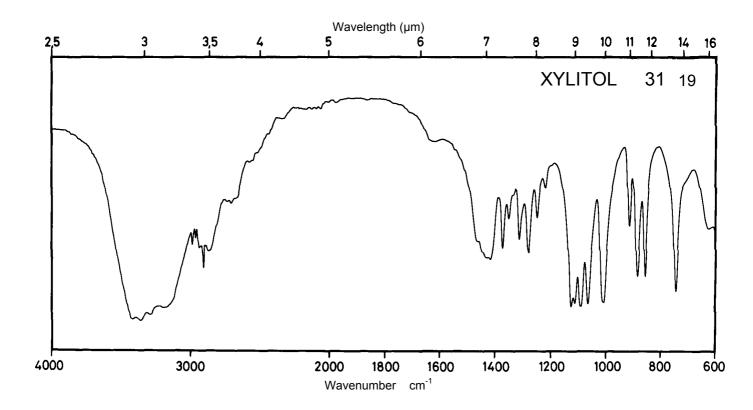


Name XYLOMETAZOLINE HYDROCHLORIDE

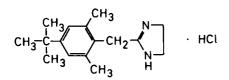


 M_{r}


280.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm		265 nm	
E ^{1%} 1cm	9.75		10.1	
ε	274		284	

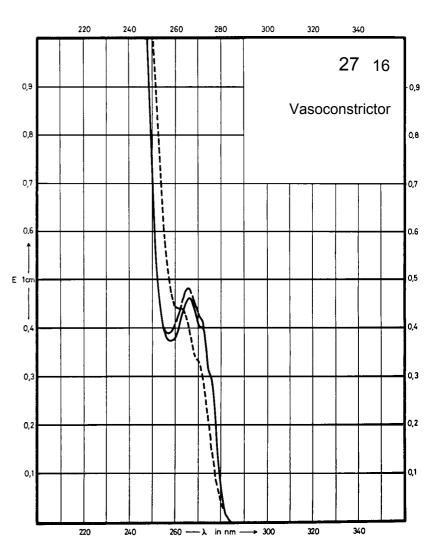

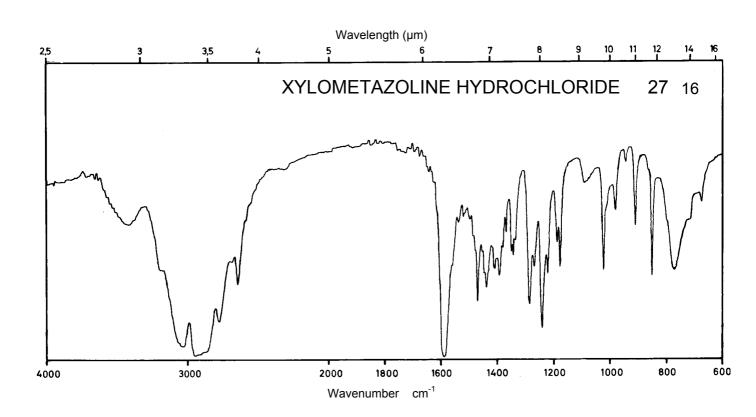





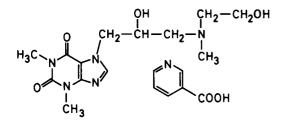

 M_{r} 434.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm		267 nm	270 nm
E ^{1%} 1cm	247		282	261
ε	10700		12300	11350

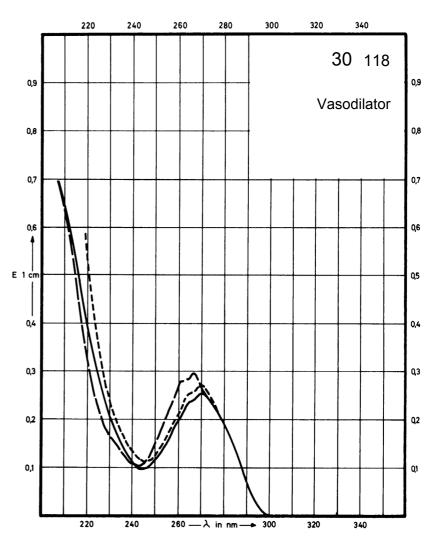

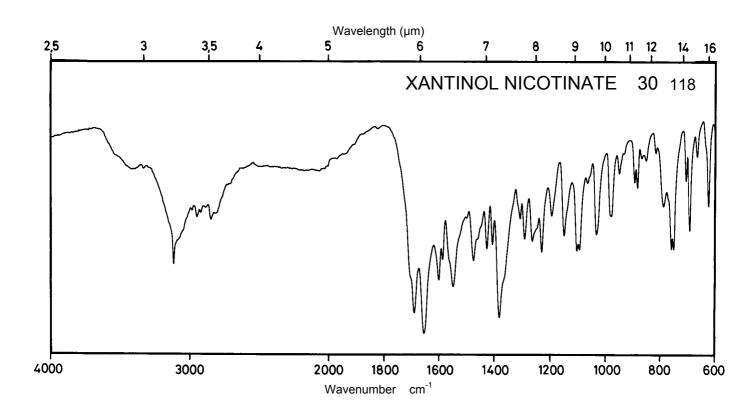



Name WARFARIN SODIUM SALT

 M_{r}


330.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm			310 nm
E ^{1%} 1cm	365			410
3	12060			13540

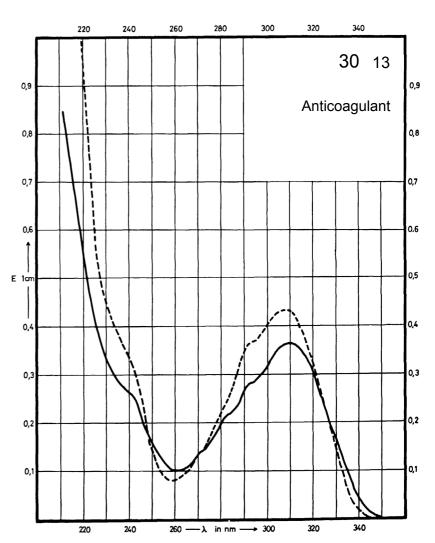

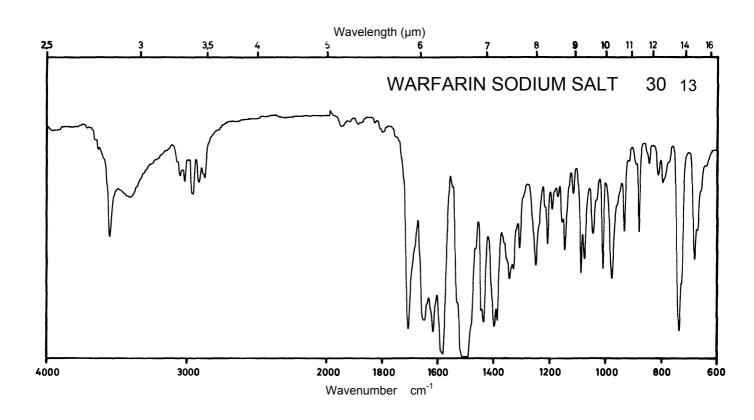




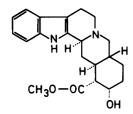


 M_r 354.5


Concentration 0.5 mg / 100 ml 2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm 226 nm		272 nm 220 nm	281 nm
E ^{1%} 1cm	219 1080		210 1110	205
ε	7780 38290		7450 39350	7260

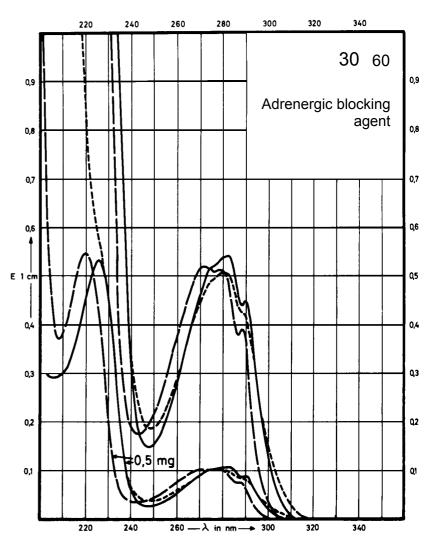

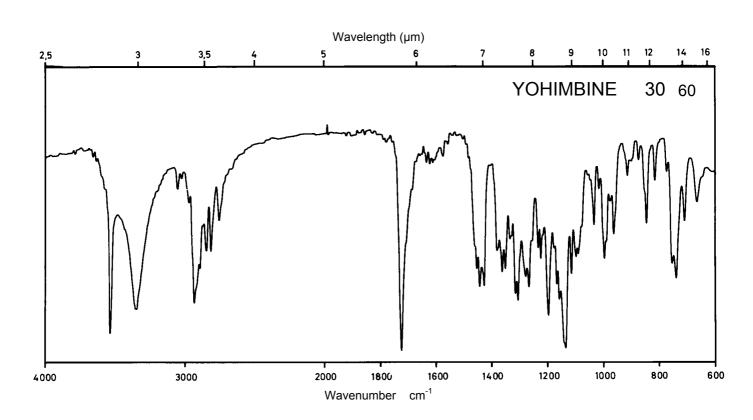


Name YOHIMBINE HYDROCHLORIDE

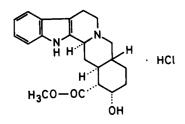


 M_r 390.9


Concentration 0.5 mg / 100 ml 2.6 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 272 nm 271 nm 279 nm absorption E^{1%} 1cm 200 187 184 ε 7800 7380 7200

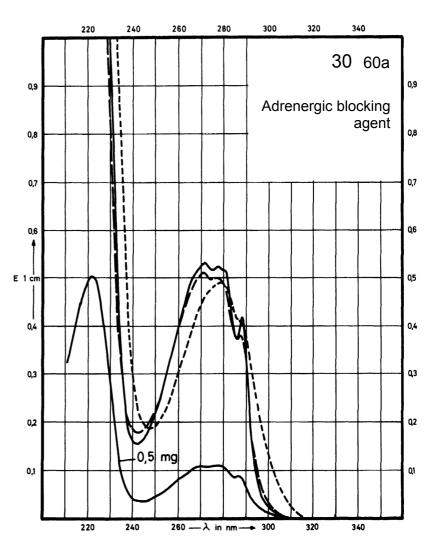

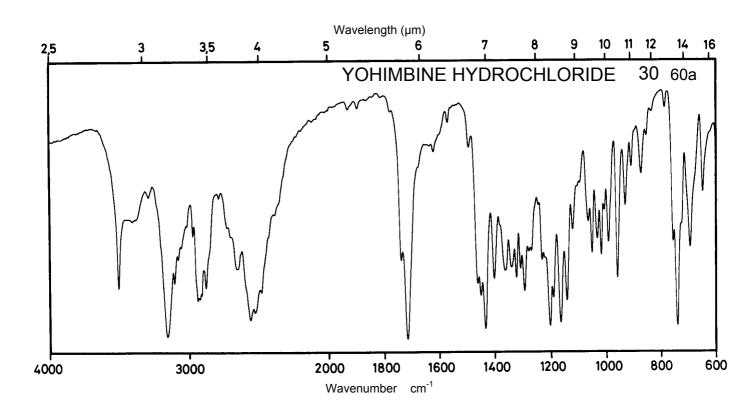




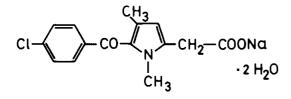


 M_r


349.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	324 nm	330 nm	324 nm	330 nm
	254 nm	260 nm	262 nm	260 nm
E ^{1%}	432	435	410	435
1cm	400	352	349	352
ε	13600	13700	12900	13700
	12600	11000	11000	11000

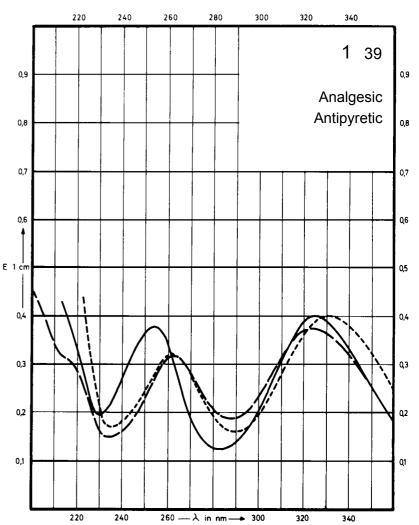

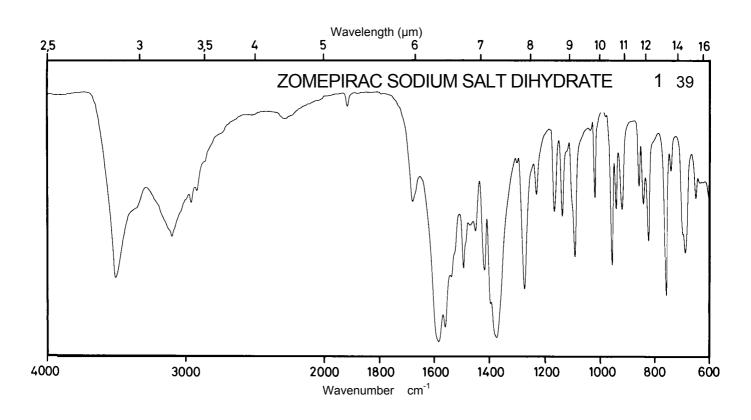



Name ZERANOL

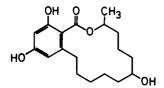


 M_r 322.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	304 nm 265 nm 218 nm	260 nm	260 nm	303 nm
E ^{1%} 1cm	164 404 767	202	205	153
ε	5300 13000 24740	6530	6600	4950

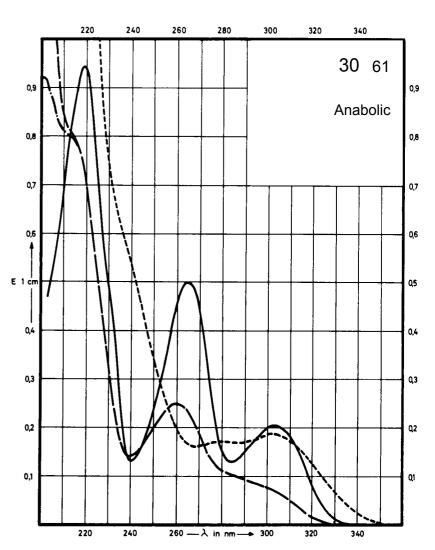

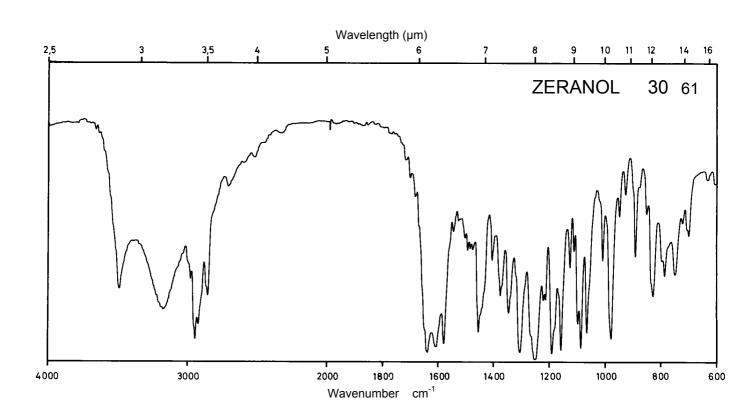




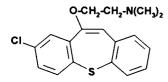


 M_{r}


331.9

[
Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
-,				
Maximum of absorption	296 nm 263 nm		295 nm 261 nm	
E ^{1%} 1cm	186 556		193 530	
ε	6200 18400		6400 17600	

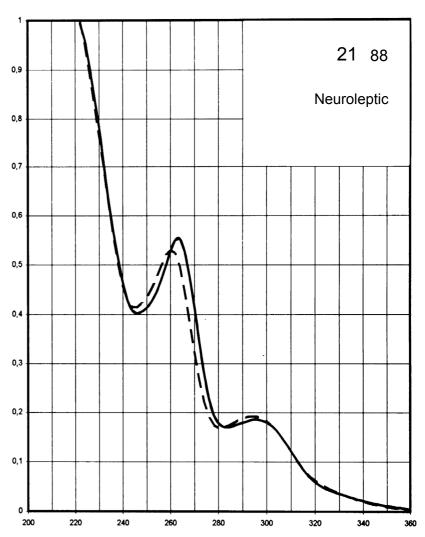

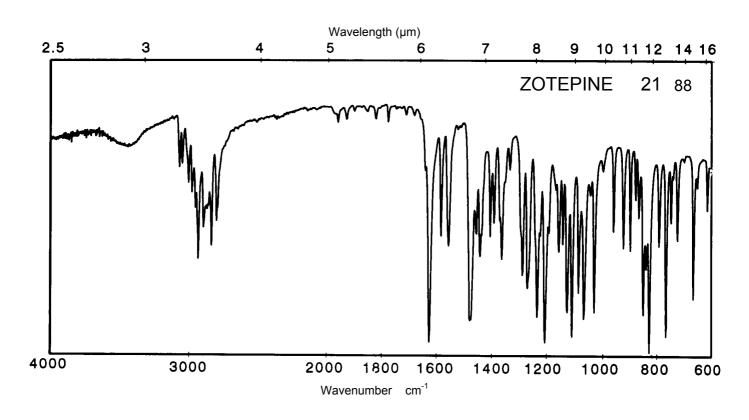


Name ZUCLOPENTHIXOL DIHYDROCHLORIDE

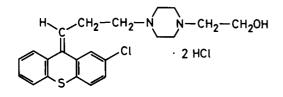


 M_{r}


473.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	326 nm 269 nm 229 nm		324 nm 267 nm 229 nm	322 nm 268 nm 227 nm
E ^{1%} 1cm	65 289 714		60 278 701	53 253 657
3	3060 13700 33800		2800 13200 33200	2500 12000 31150

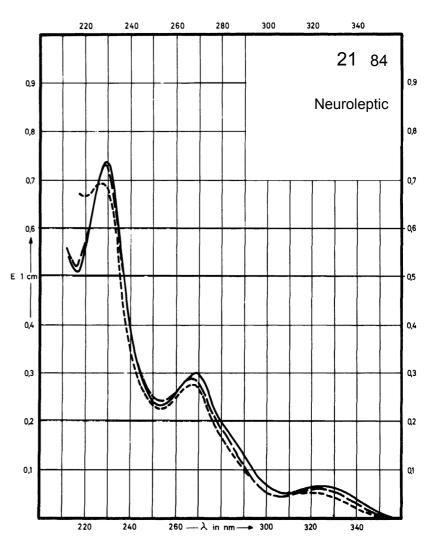

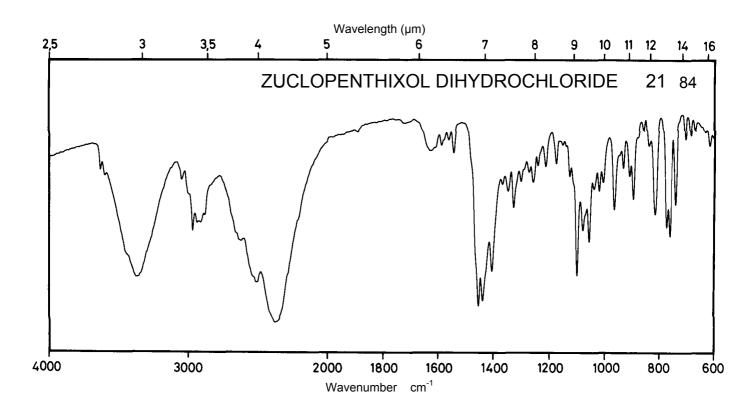

Name ZIMELIDINE DIHYDROCHLORIDE

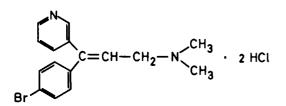



 M_r

Concentration 1 mg / 100 ml


390.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm	261 nm	250 nm	258 nm
E ^{1%} 1cm	485	485	476	486
ε	18900	18900	18600	19000

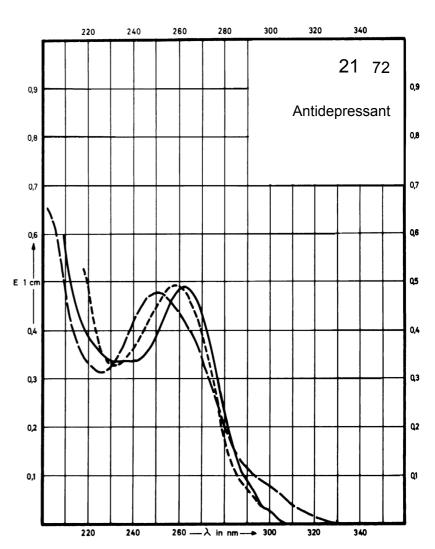

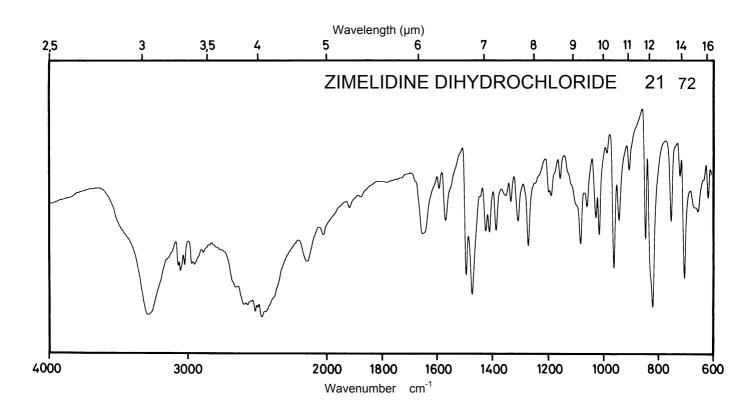




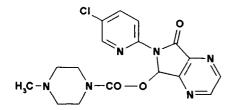


 $\,M_r\,$


388.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	303 nm	304 nm	277 nm 237 nm
E ^{1%} 1cm		362	364	199 390
3		10500	10500	5800 11300

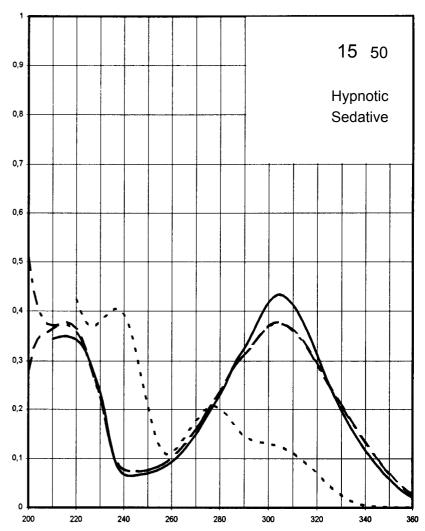

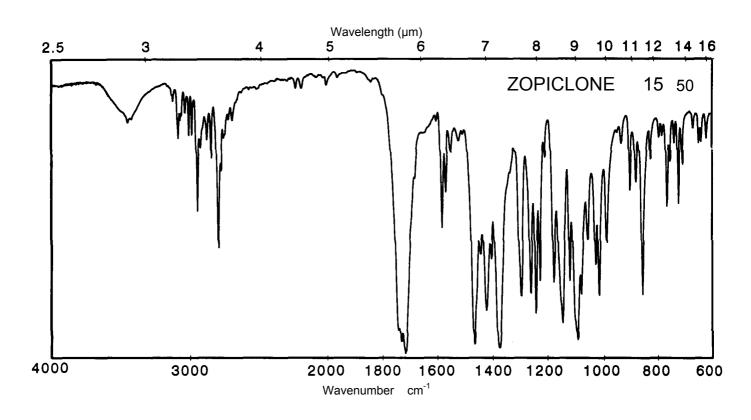


Name ZOLPIDEM TARTRATE

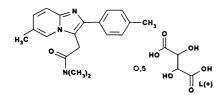


 M_{r}


382.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm 263 nm		294 nm 260 nm	
E ^{1%} 1cm	186 558		192 531	
3	7100 21300		7400 20300	

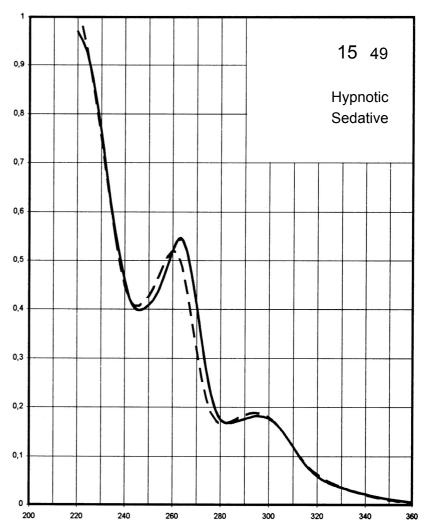

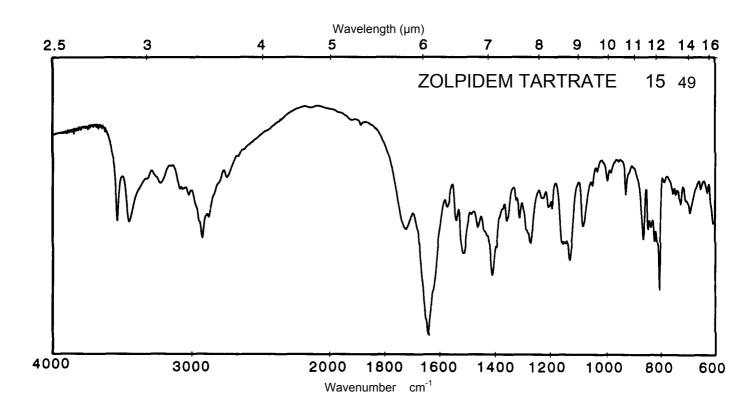



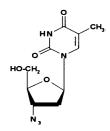


 M_{r}


267.2

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm	267 nm	266 nm	266 nm
E ^{1%} 1cm	387	381	362	379
ε	10300	10200	10200	10100

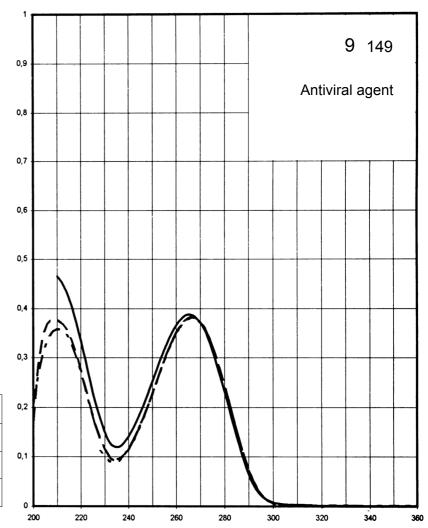

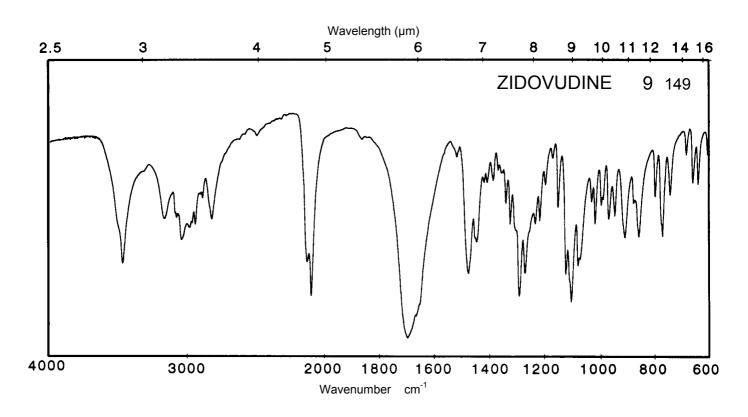




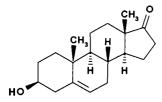


 M_{r}


288.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	294 nm			
E ^{1%} 1cm	1.6			
3	47			

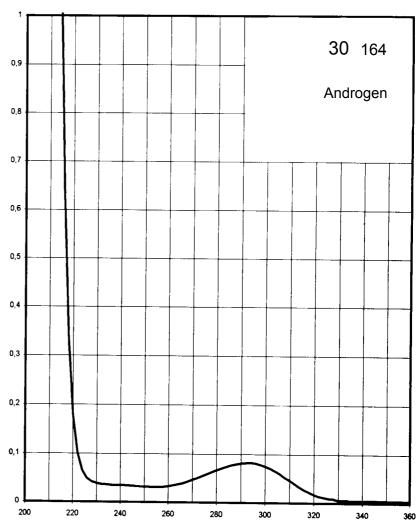

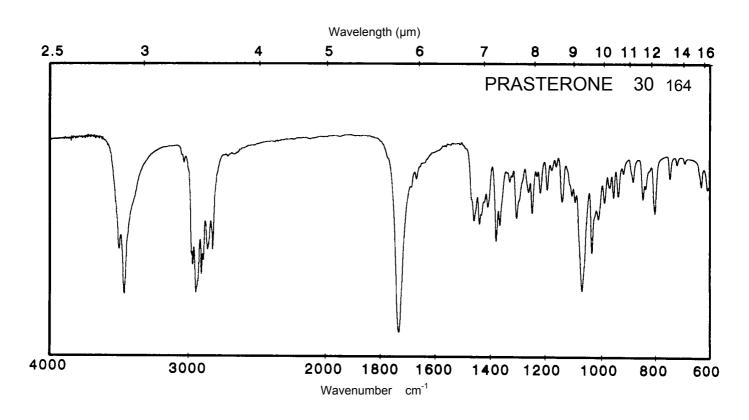



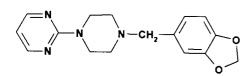



 M_r

Concentration 1 mg / 100 ml


298.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 243 nm	286 nm 242 nm	312 nm 287 nm 236 nm	286 nm 243 nm
E ^{1%} 1cm	196 939	179 823	125 204 893	177 865
3	5800 27000	5300 24500	3700 6100	5300 25800

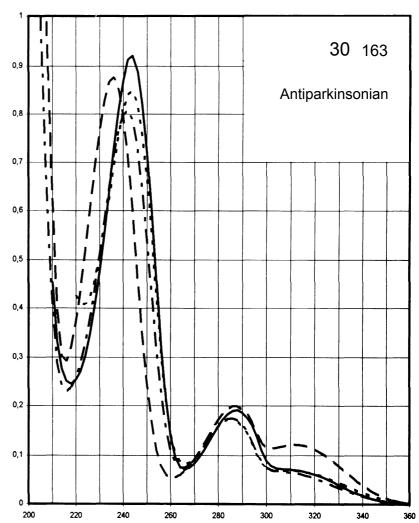

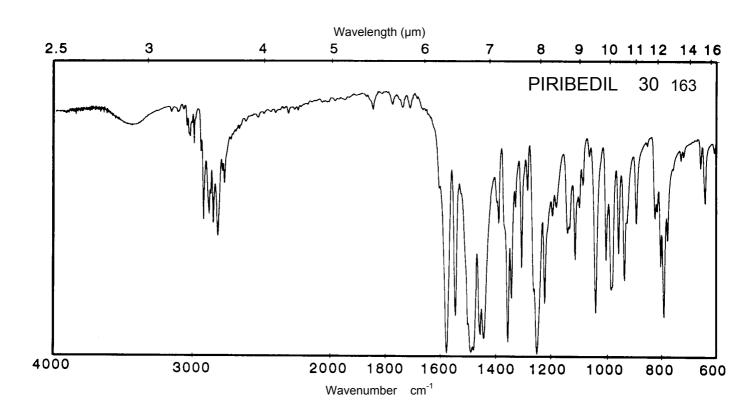





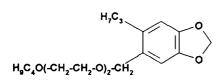


 M_{r} 338.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	290 nm	289 nm	289 nm	289 nm
	238 nm	238 nm	238 nm	236 nm
E ^{1%}	129	112	107	100
1cm	151	140	136	130
ε	4350	3800	3600	3400
	4400	4700	4600	4400

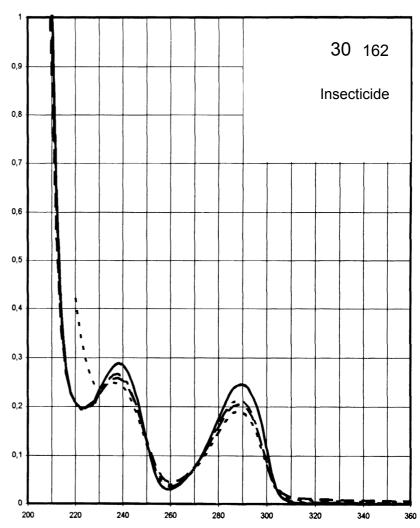

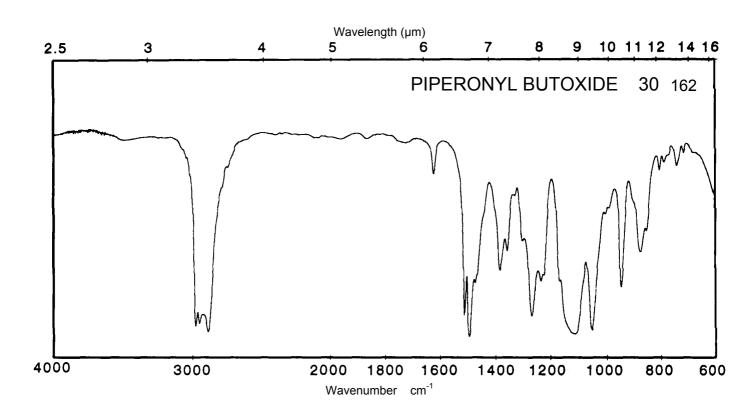




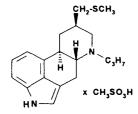


 M_{r}


410.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm 280 nm 222 nm	279 nm 222 nm	279 nm 221 nm	
E ^{1%} 1cm	153 183 865	150 725	160 804	
ε	6300 7500 35500	6150 29800	6600 33000	

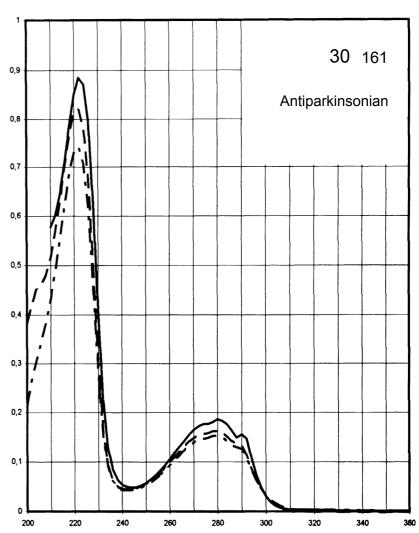

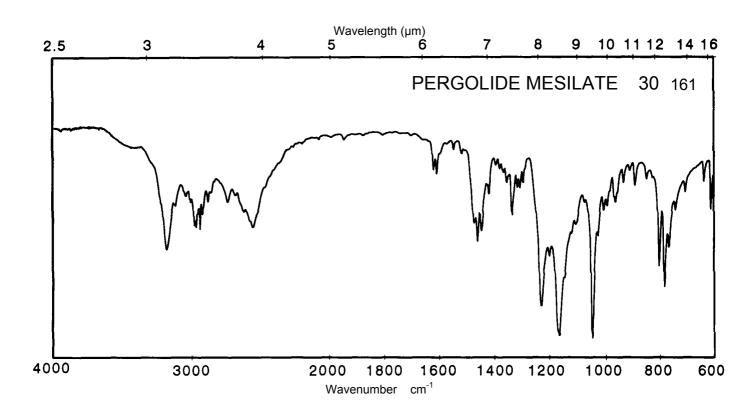





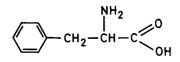

 M_r 165.2

Concentration 75 mg / 100 ml

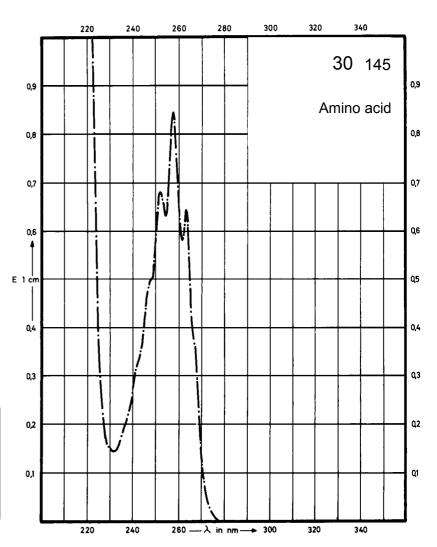

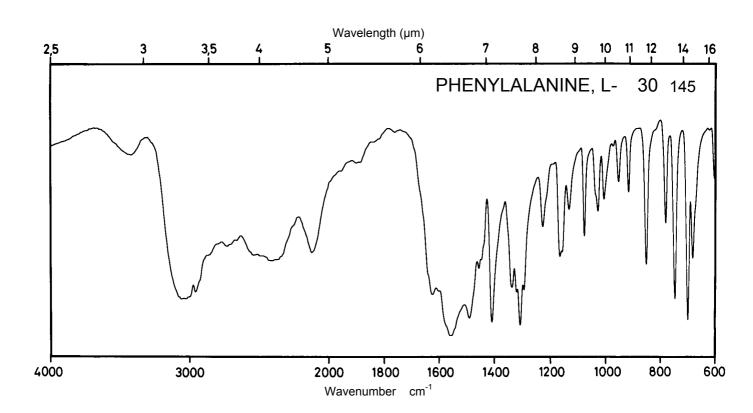








 M_r 115.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

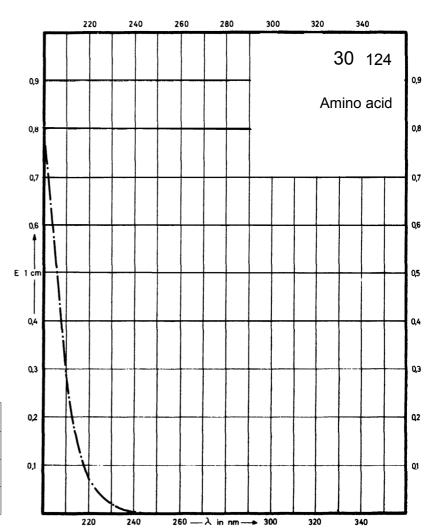

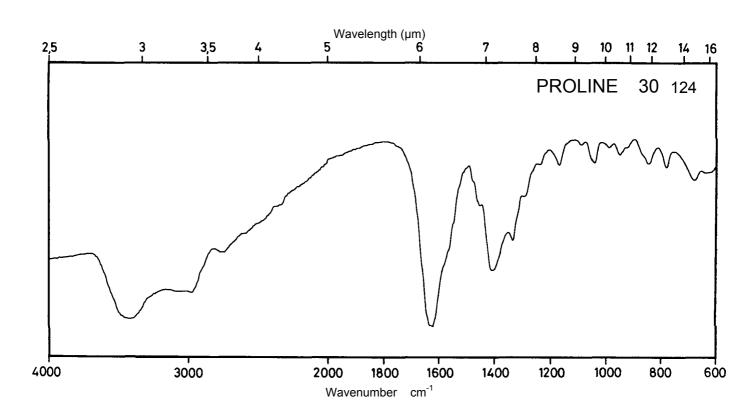




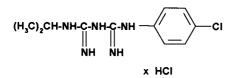


 M_r


290.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm	254 nm 232 nm		239 nm
E ^{1%} 1cm	715	536 509		498
3	20700	15600 14800		14400

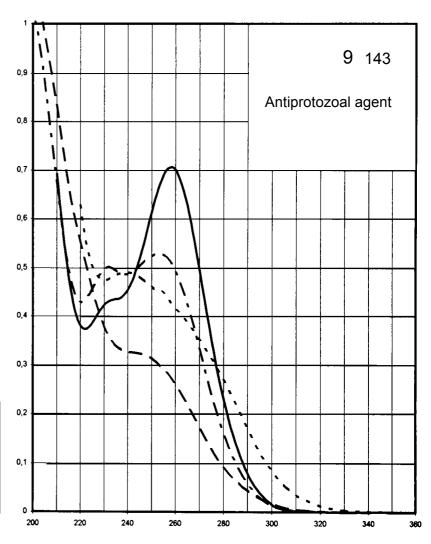

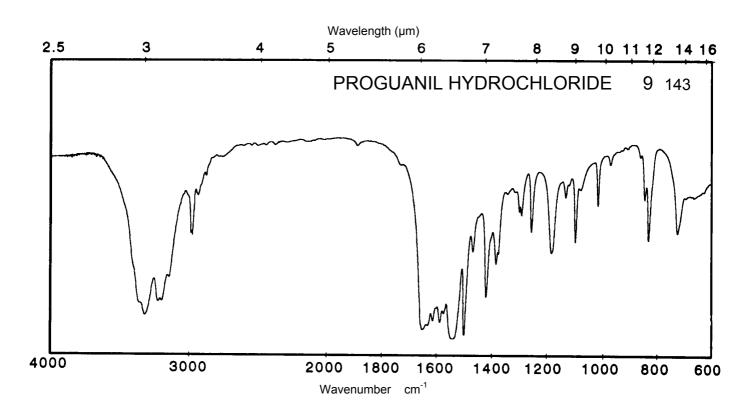





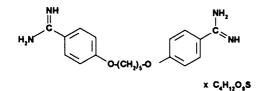

592.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm	262 nm	262 nm	248 nm
E ^{1%} 1cm	509	463	464	431
3	30200	27400	27500	25500

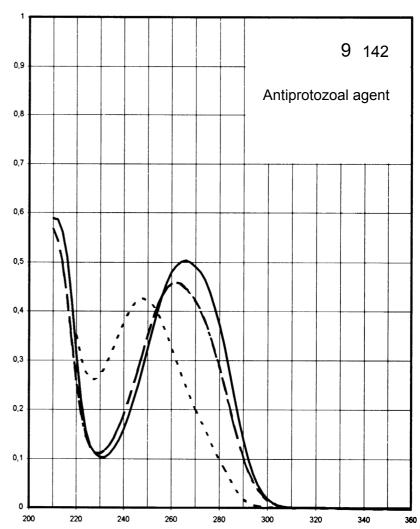


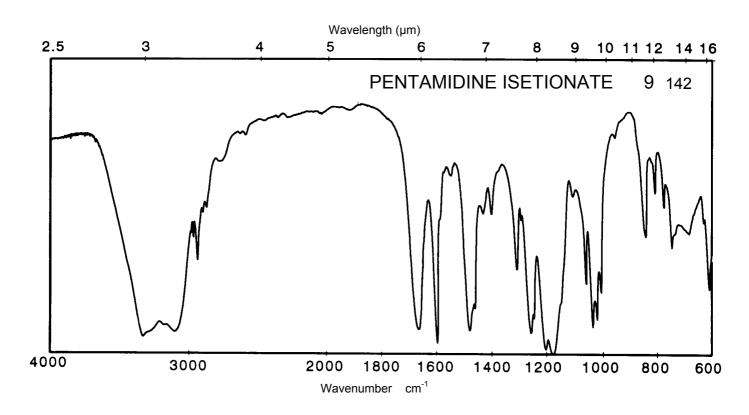


Name PEFLOXACIN MESILATE


 M_{r} 465.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	317 nm	315 nm	316 nm	323 nm
	278 nm	276 nm	277 nm	272 nm
E ^{1%}	278	266	266	306
1cm	955	911	950	814
3	12900	12400	12400	14200
	44400	42400	44200	37900

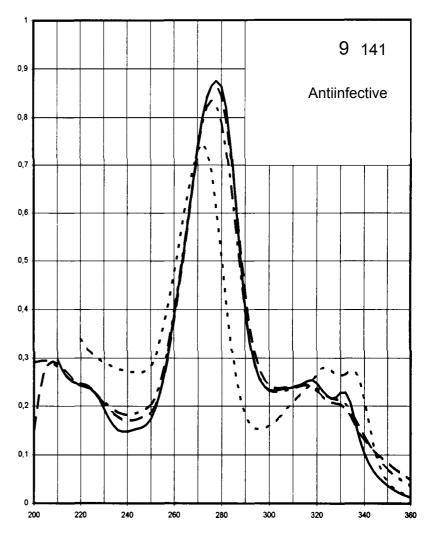

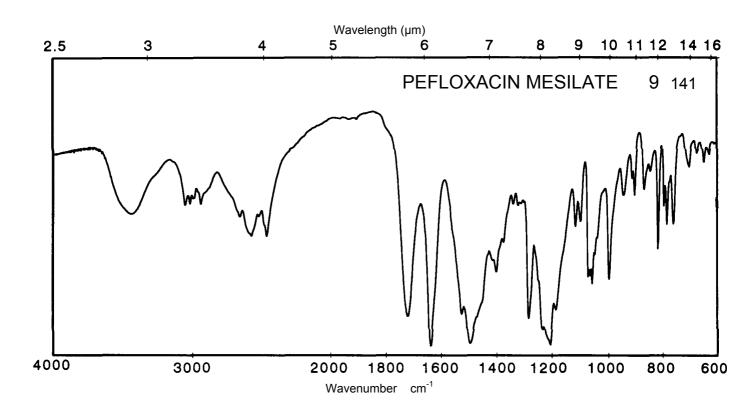




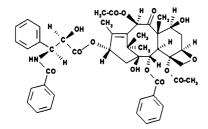



 M_{r}


853.9


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	228 nm		231 nm	
E ^{1%} 1cm	403		333	
3	34400		28400	

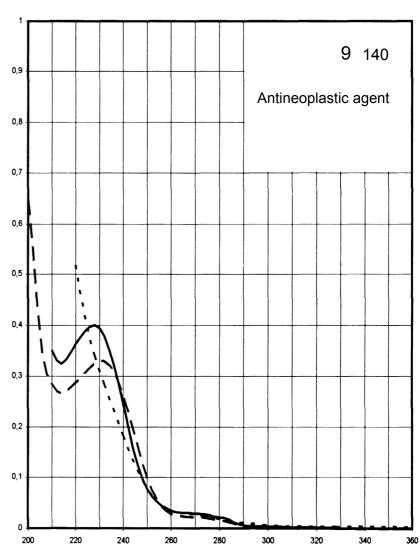

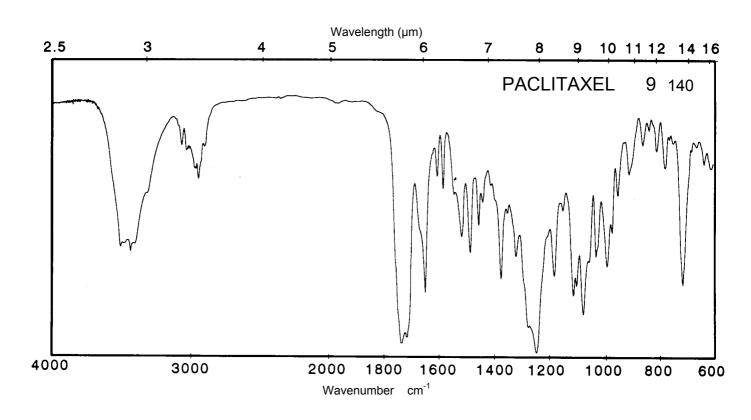

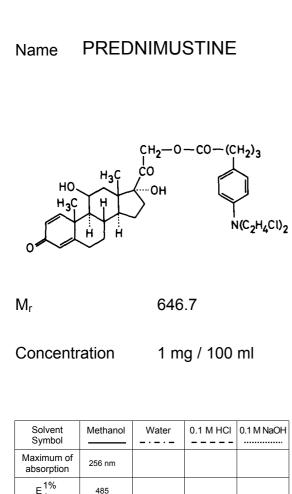

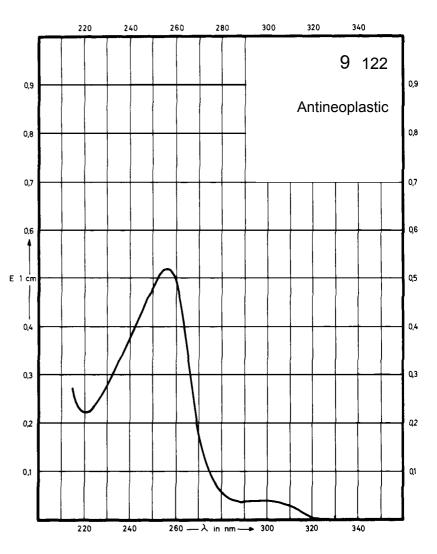


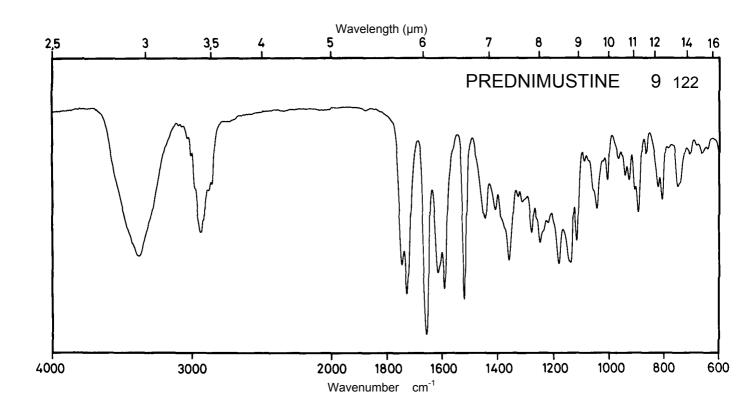






 $HO-CH_2-(CH_2-O-CH_2)_n-CH_2-OH$

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

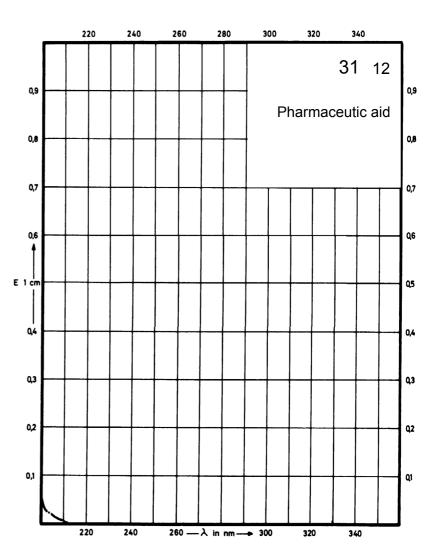


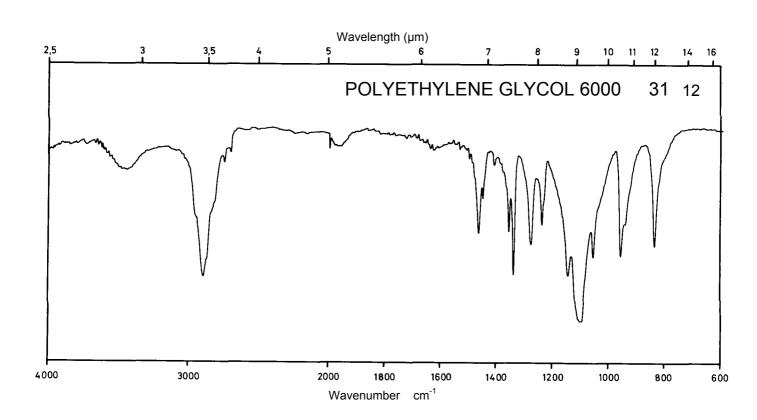


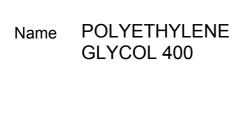

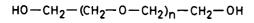

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

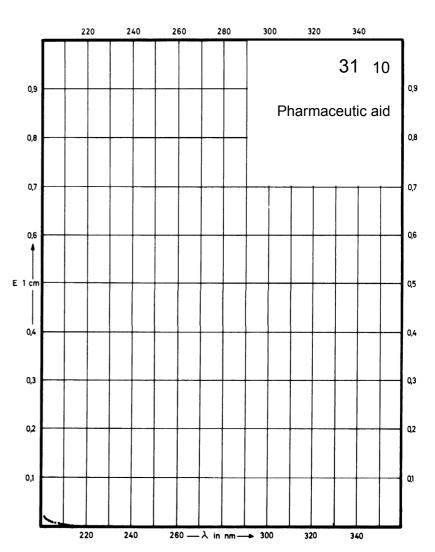

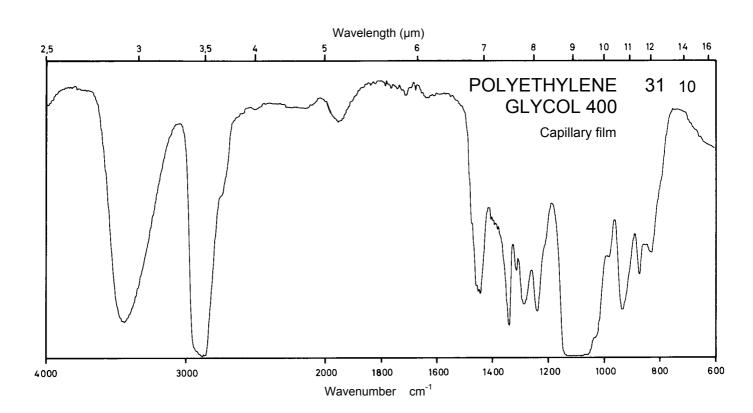





СН₃—СН—СН₂—ОН ОН


 M_r 76.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

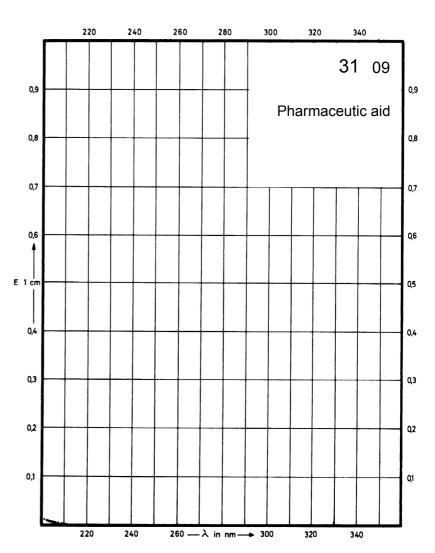


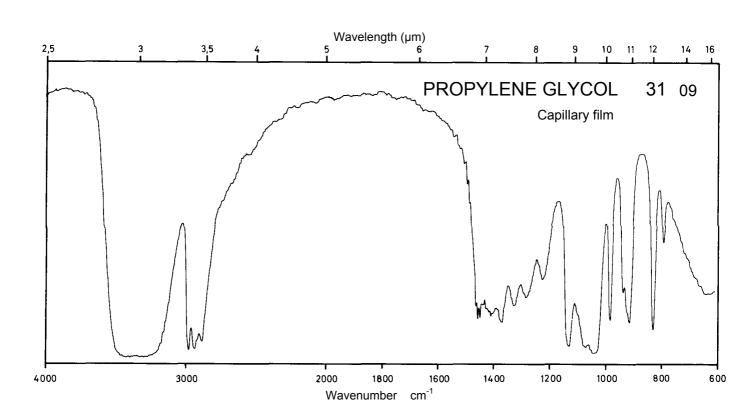





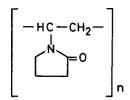


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
ε				

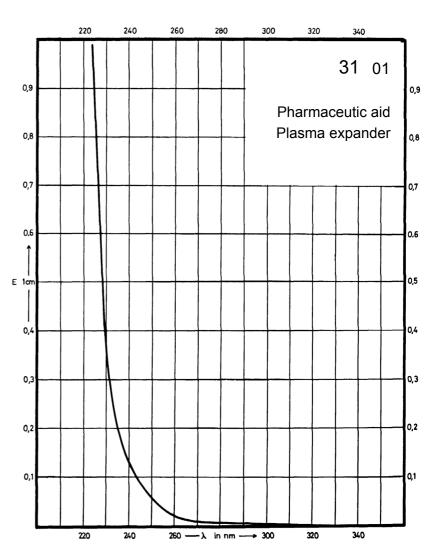


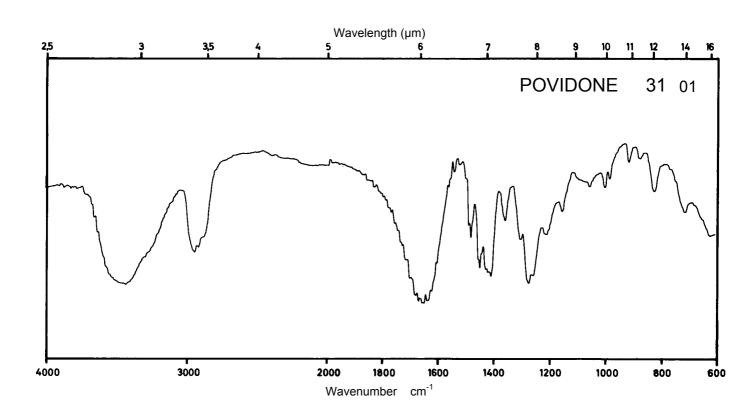

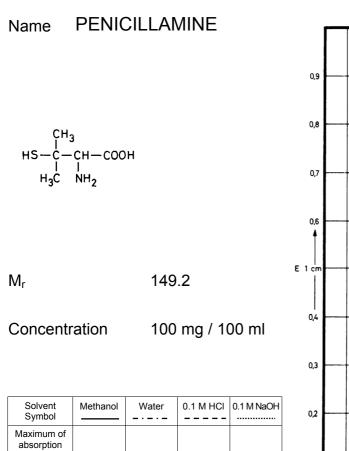

ε

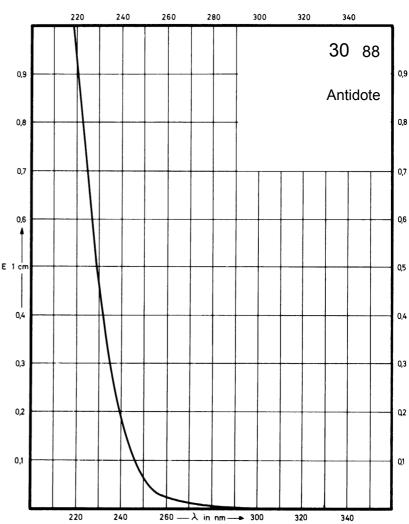

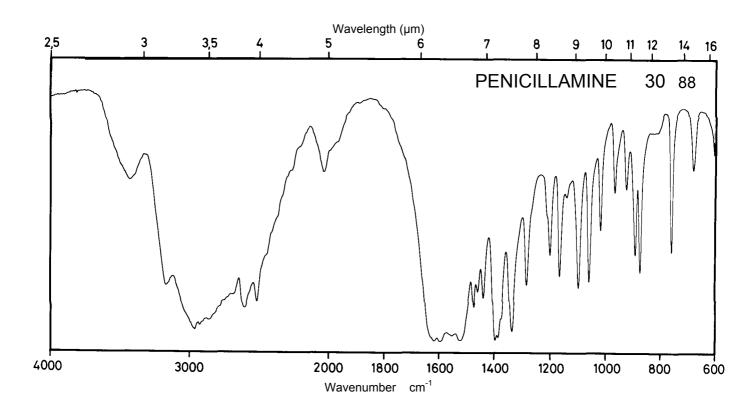








109.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm	261 nm	261 nm	261 nm
E ^{1%} 1cm	256	268	494	262
ε	2790	2920	5390	2860

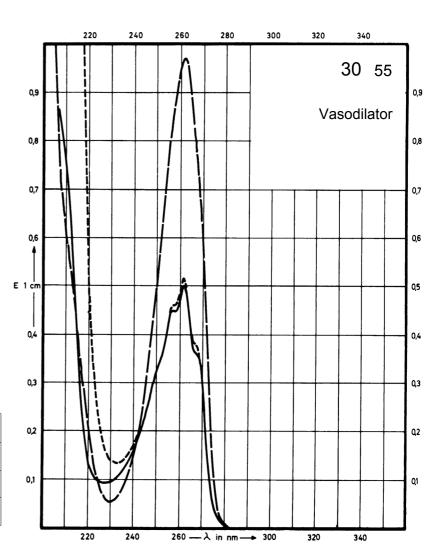

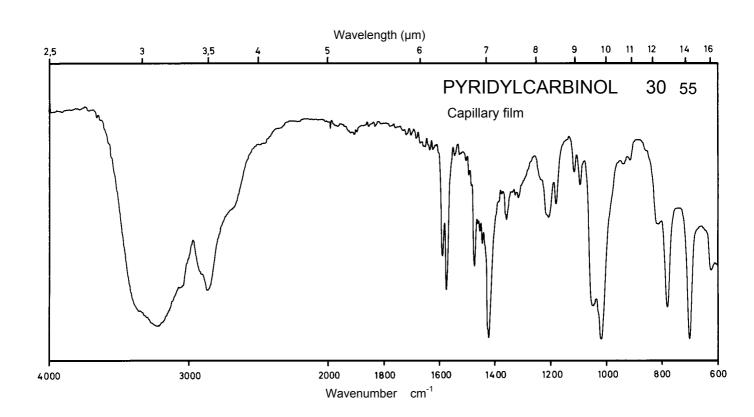





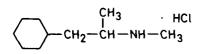

191.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
ε				

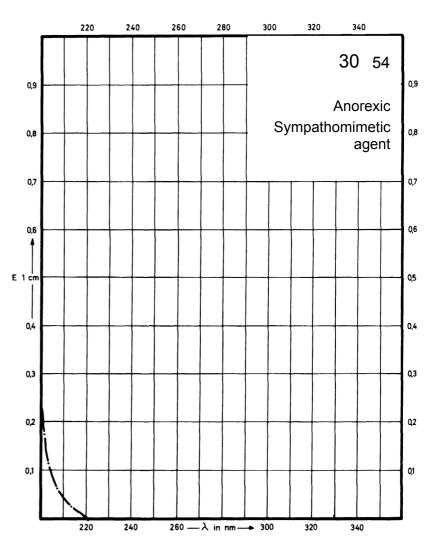

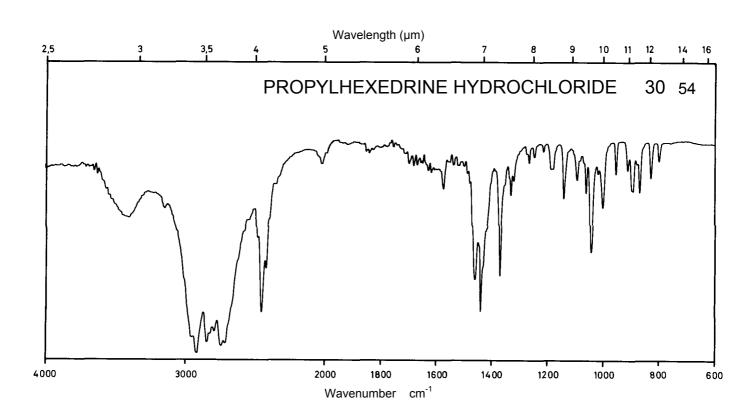



Name PHENGLUTARIMIDE HYDROCHLORIDE

 M_r


324.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm	264 nm 258 nm 252 nm		
E ^{1%} 1cm	8.5 10.6 9.7	8.9 11.8 11.5		
ε	275 344 315	290 383 374		

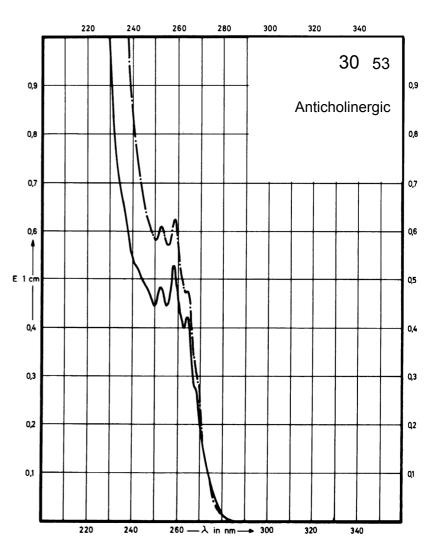

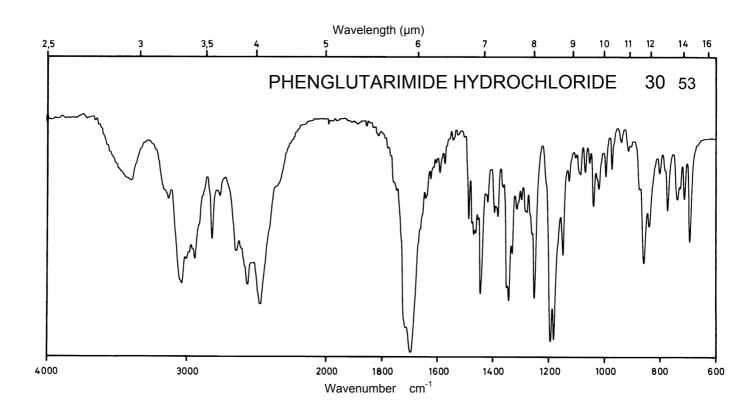





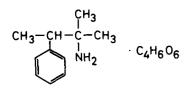


 M_r 313.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm		264 nm 258 nm 252 nm	
E ^{1%} 1cm	4.26 5.51 4.29		4.23 5.65 4.50	
ε	134 173 134		133 177 141	

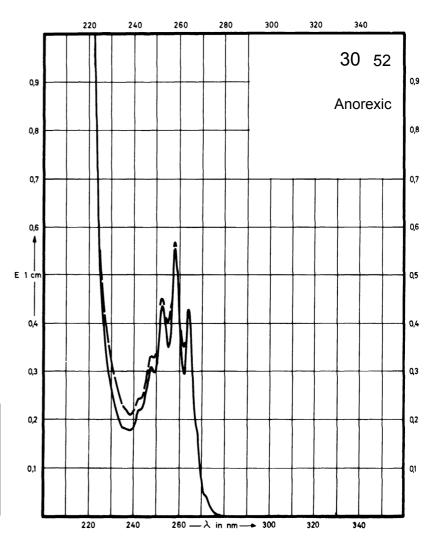

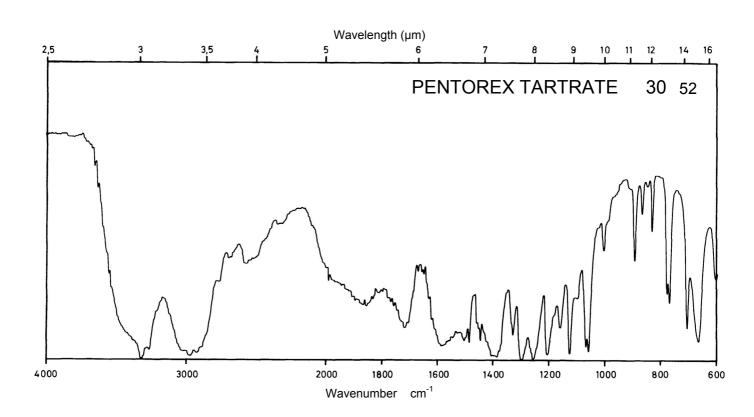





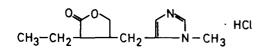



244.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	216 nm		215 nm	
E ^{1%} 1cm	243		241	
ε	5940		5900	

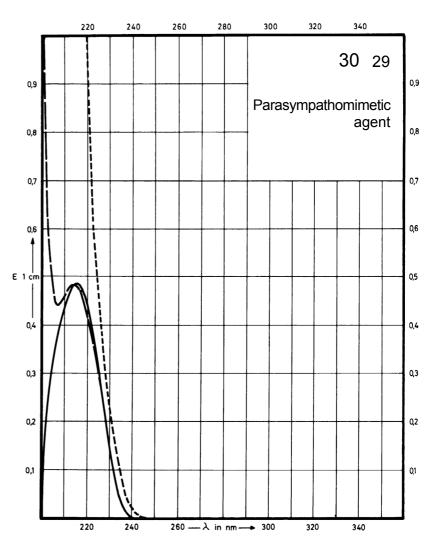


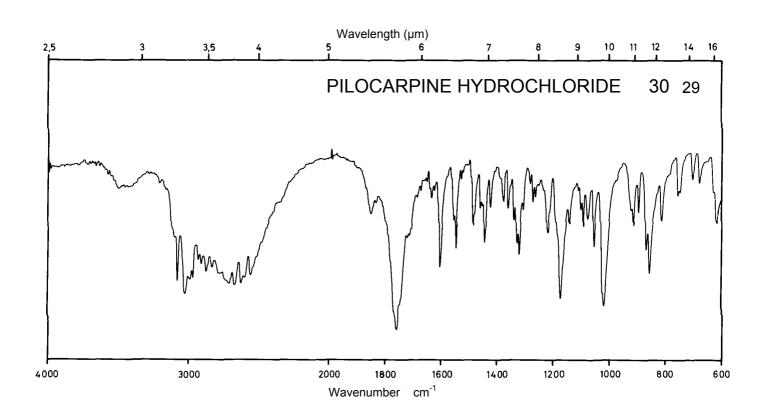

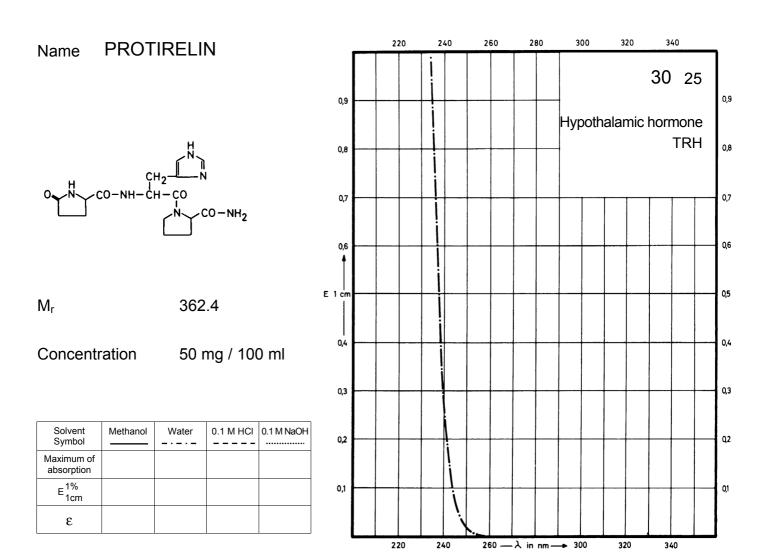

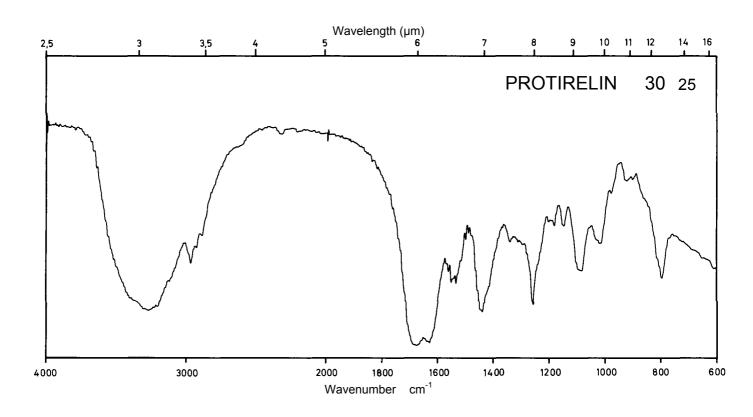





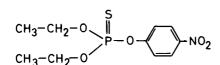


291.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm			
E ^{1%} 1cm	335			
ε	9760			

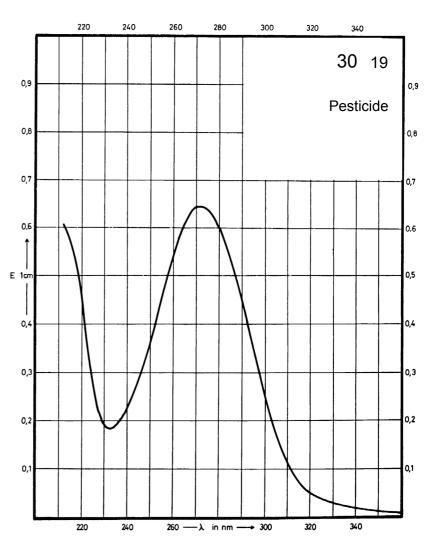

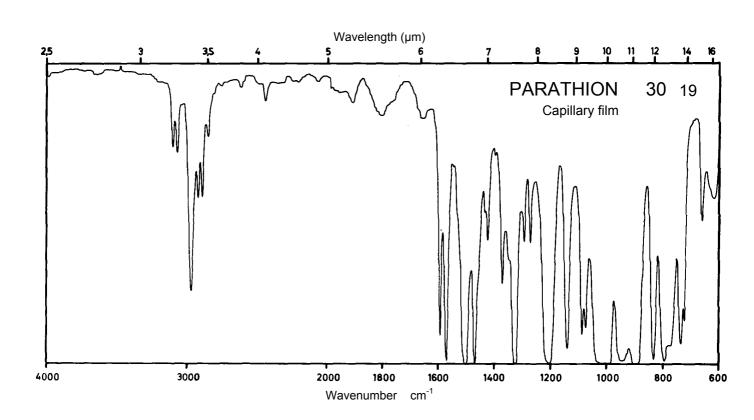




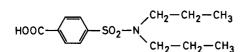


 M_r


285.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	248 nm			243 nm
E ^{1%} 1cm	332			348
ε	9470			9930

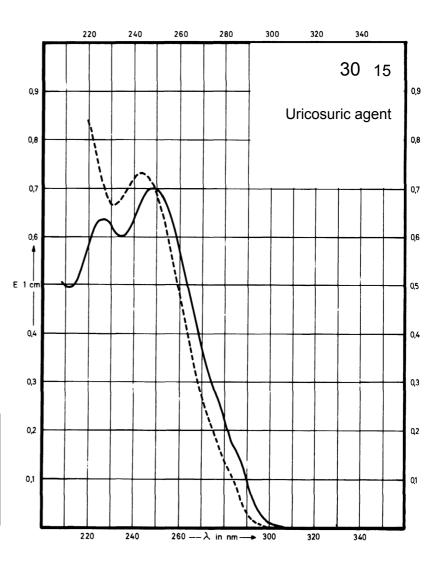

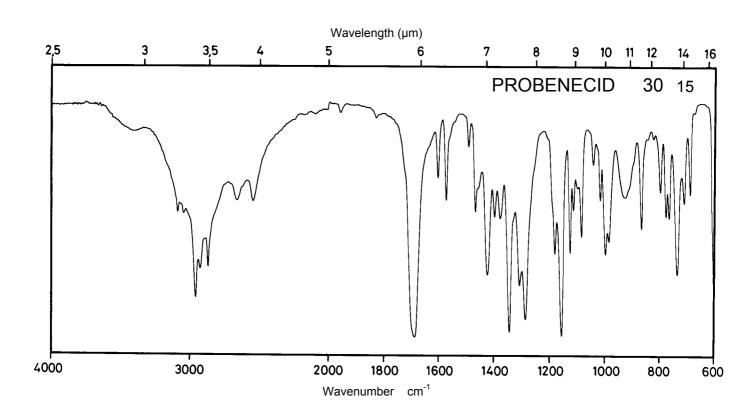



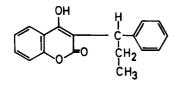



280.3

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm 285 nm			310 nm
E ^{1%} 1cm	440 407			540
3	12330 11410			15140

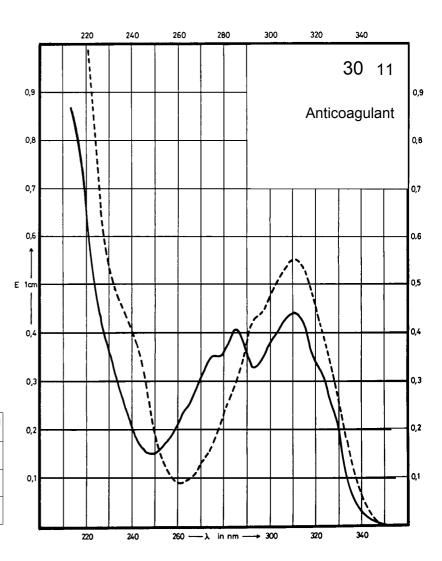

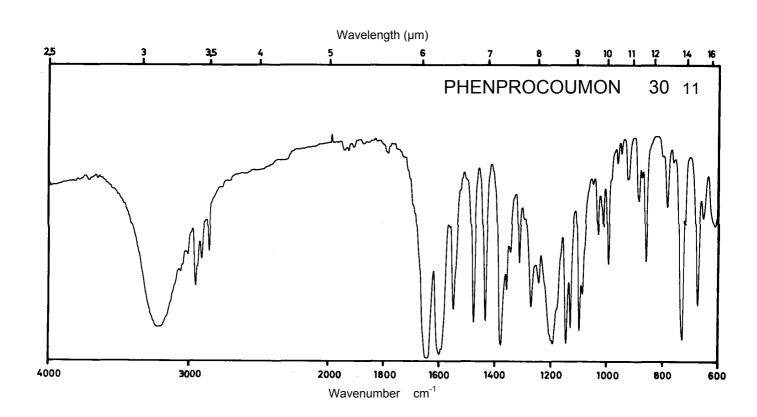





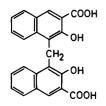


 M_r 388.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	369 nm 288 nm			363 nm 287 nm
E ^{1%} 1cm	173 266			171 252
ε	6720 10330			6640 9790

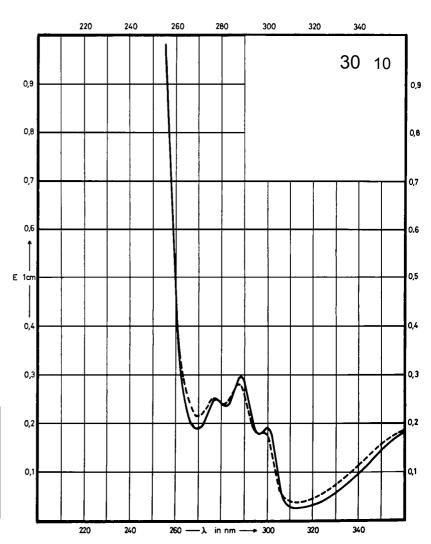

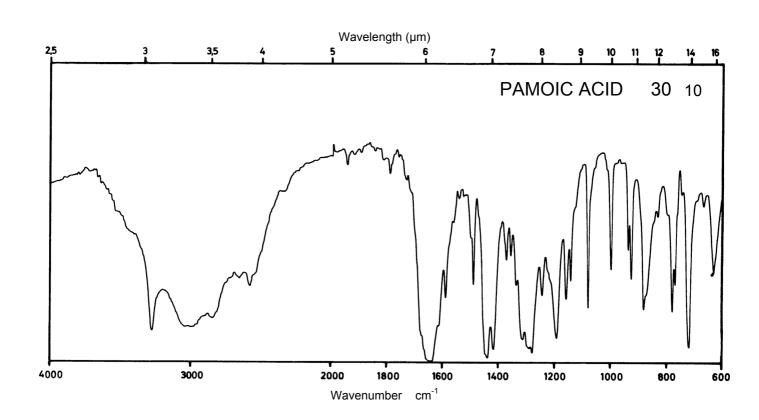





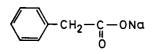


158.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm		256 nm	257 nm
E ^{1%} 1cm	11.8		10.8	12.3
ε	190		170	190

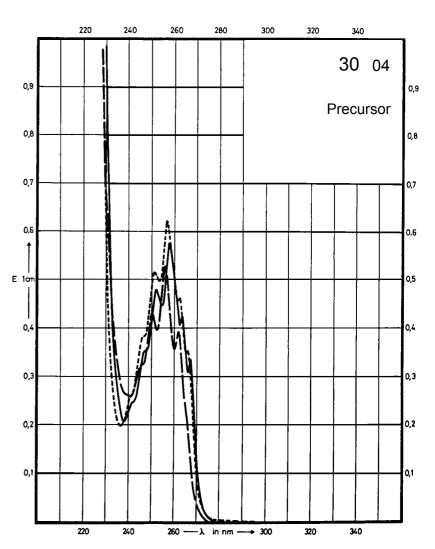

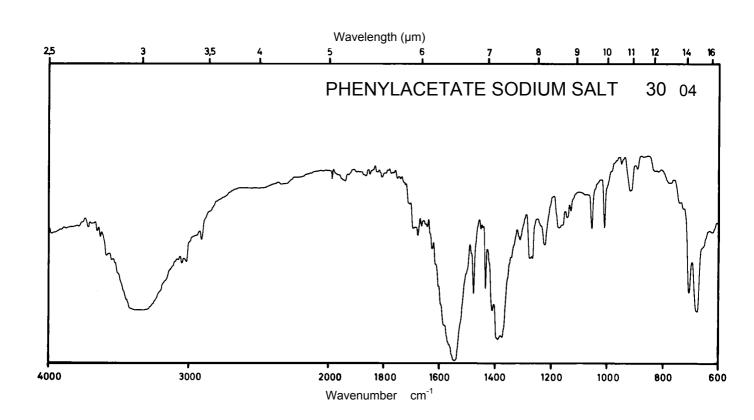




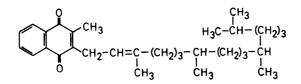


 M_r 450.7

Concentration 1.7 mg / 100 ml 8.3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	330 nm 271 nm 249 nm			
E ^{1%} 1cm	67 348 388			
ε	3030 15700			

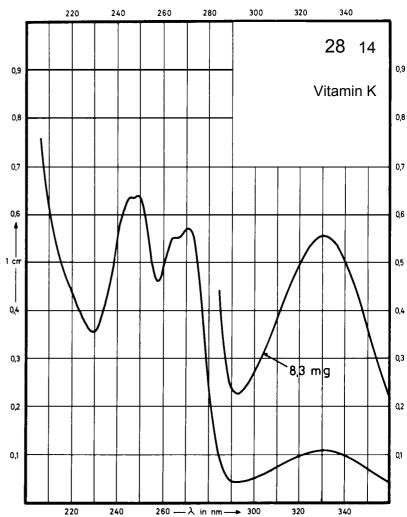

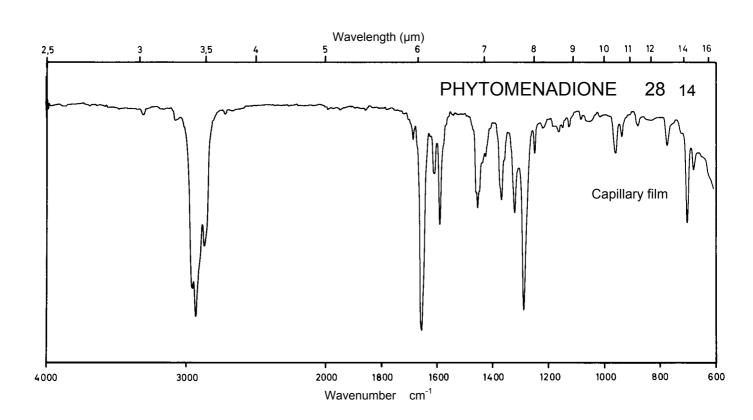



Name PYRIDOXINE HYDROCHLORIDE

 M_{r}

205.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	292 nm		290 nm	309 nm 244 nm
E ^{1%} 1cm	412		430	344 313
ε	8470		8840	7070 6440

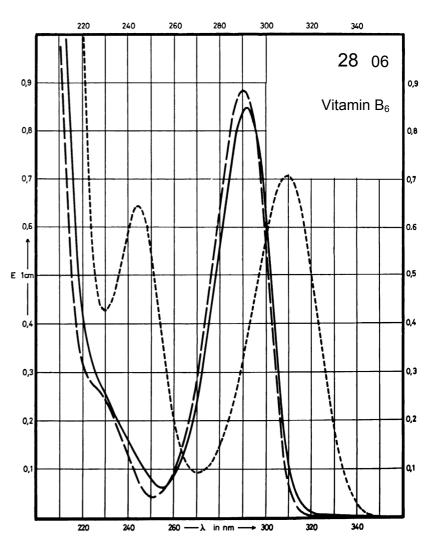

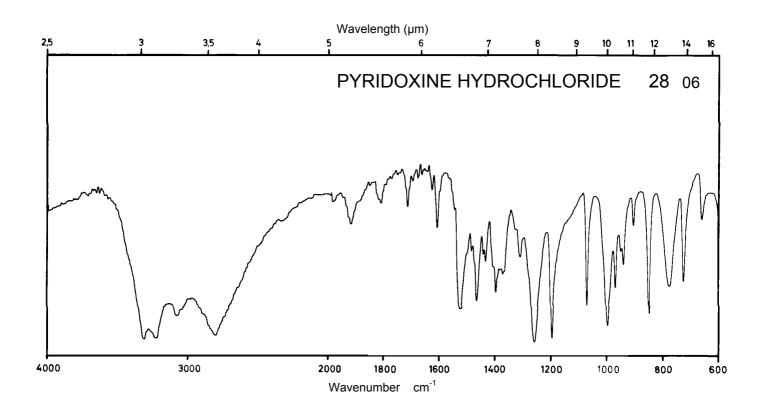



Name PHENYLPROPA-NOLAMINE HYDROCHLORIDE

 M_r


187.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 257 nm 251 nm	262 nm 256 nm 250 nm		
E ^{1%} 1cm	7.3 9.4 7.4	7.6 9.9 7.7		
ε	136 177 139	143 186 145		

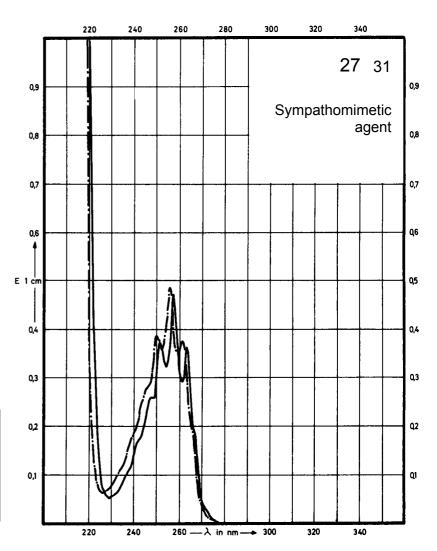

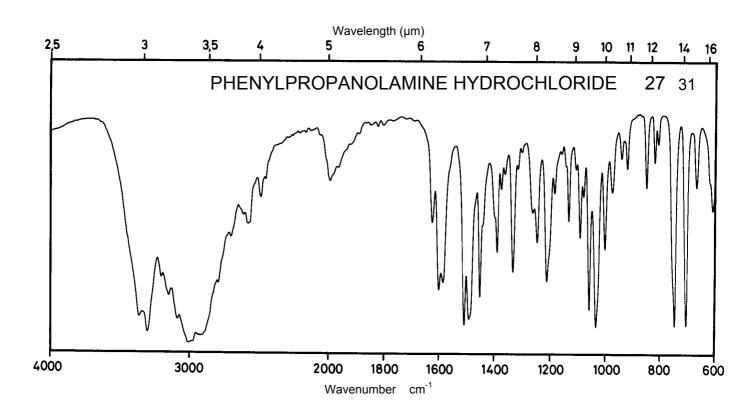





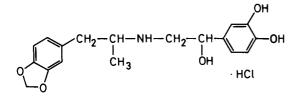

 M_r 367.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm 231 nm		282 nm 225 nm	Decom- position observed
E ^{1%} _{1cm}	199 254		167 242	
3	7300 9350		6150 8900	

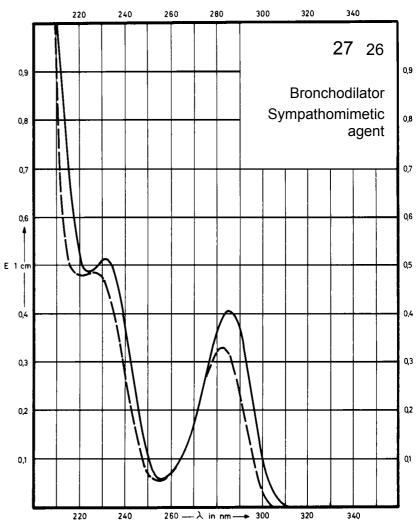

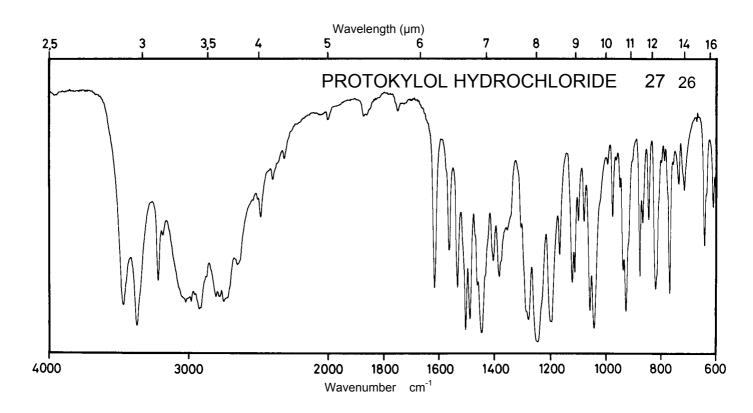





$$\left[\begin{array}{c} \text{HO-} \\ \text{CH}_2 - \text{CH-} \\ \text{CH}_3 \end{array} \right]_2 \text{SO}_4^{2^-}$$

428.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm 224 nm		275 nm 222 nm	295 nm 238 nm
E ^{1%} 1cm	78 397		65 345	116 494
ε	3340 17020		2760 14790	4970 21170

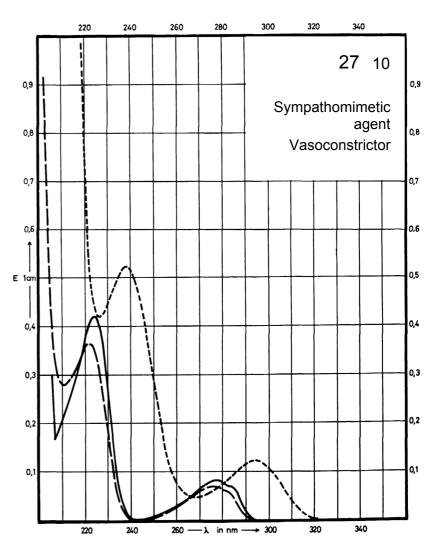

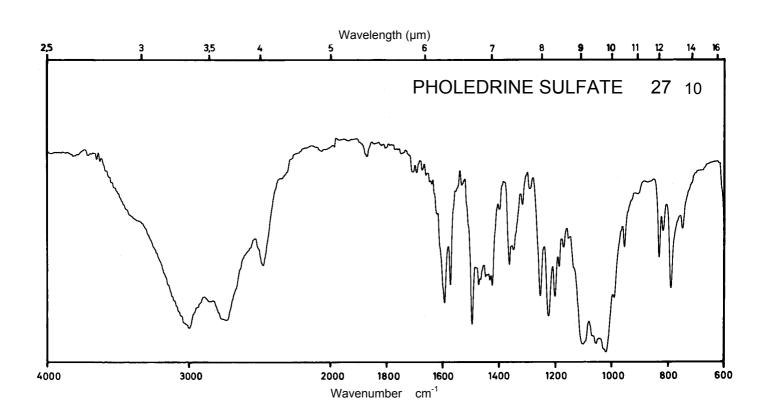



Name PHENYLEPHRINE HYDROCHLORIDE

 M_r 203.7


Concentration 1.1 mg / 100 ml 2.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm		272 nm	290 nm 236 nm
E ^{1%} 1cm	106		92	152 446
ε	2160		1870	3100 9090

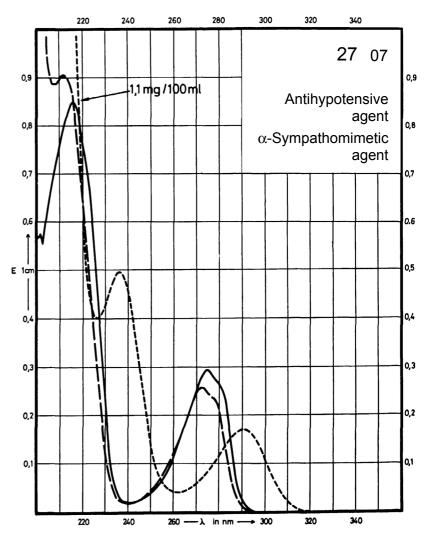

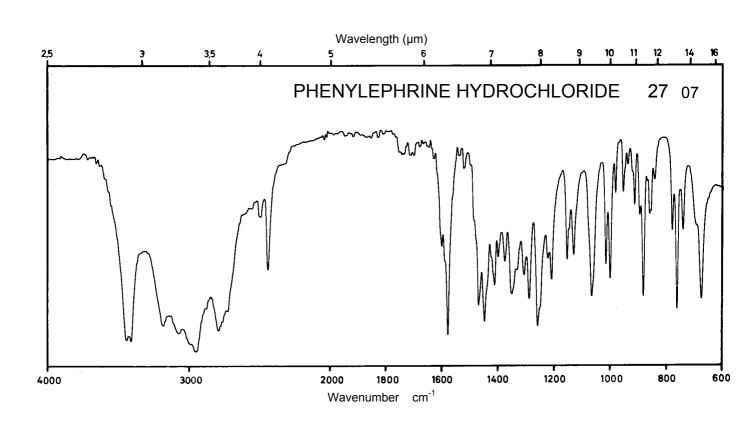




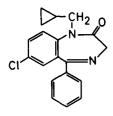


 M_r 324.8

Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 357 nm 286 nm 240 nm 118 429 Maximum of 313 nm 228 nm 310 nm absorption E^{1%} 1cm 61 920 61 820 3830 1980 29880 13930 26630 ε 1980

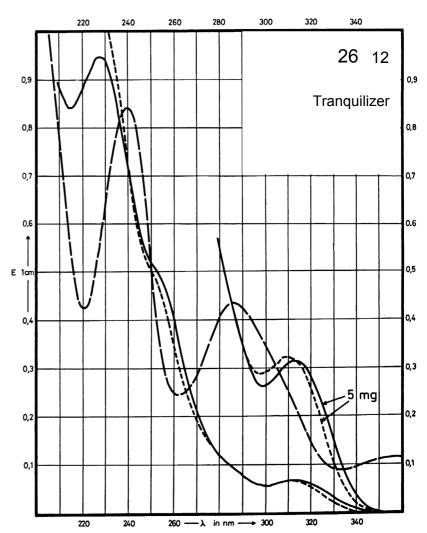

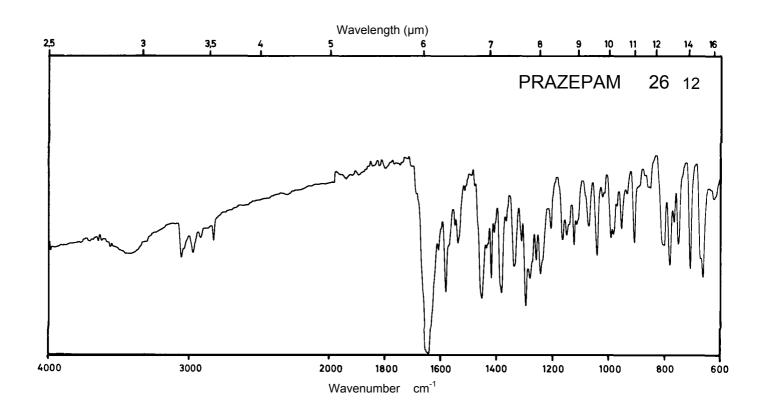

Name PYRITHYLDIONE

0 C₂H₅


 M_{r}


167.2

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	302 nm		306 nm	365 nm
E ^{1%} 1cm	500		493	642
ε	8360		8240	10730

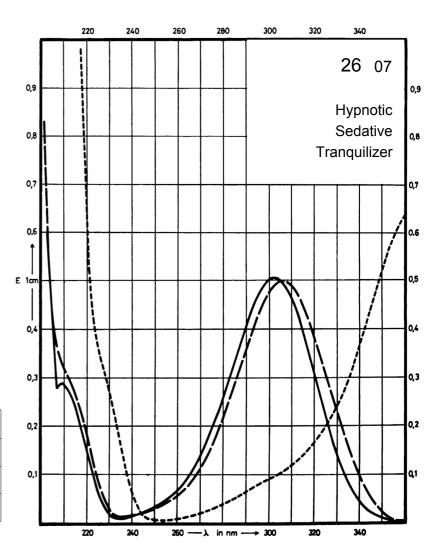

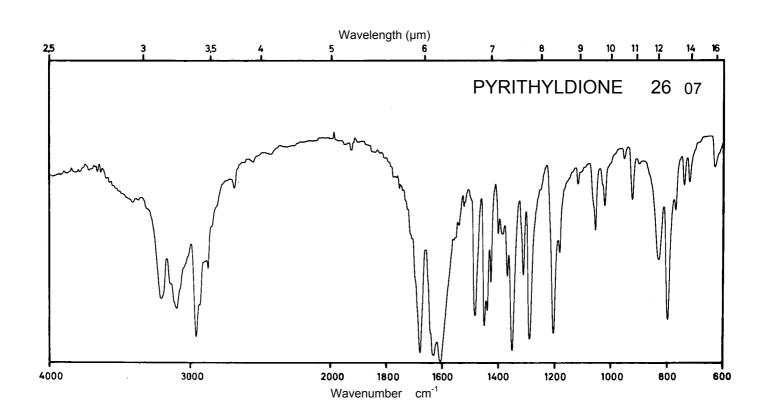




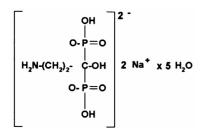


 M_r

369.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

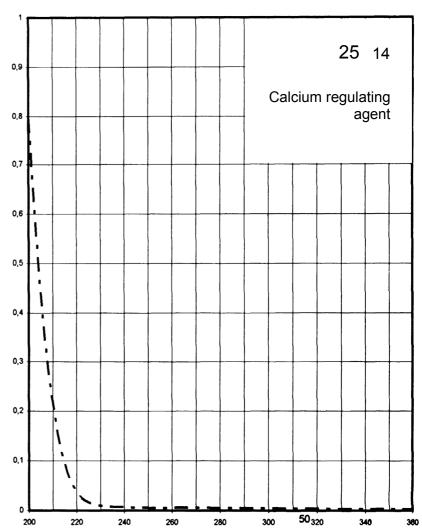

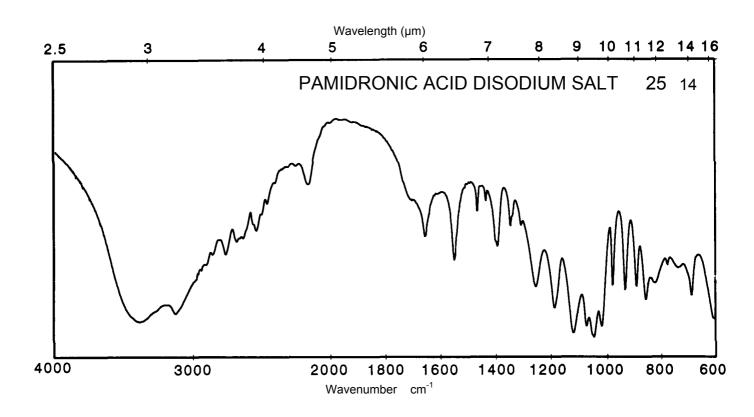



Name PROPYLTHIOURACIL

 M_{r}


170.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 215 nm		275 nm 215 nm	260 nm
E ^{1%} 1cm	957 952		987 975	704
ε	16300 16200		16800 16600	12000

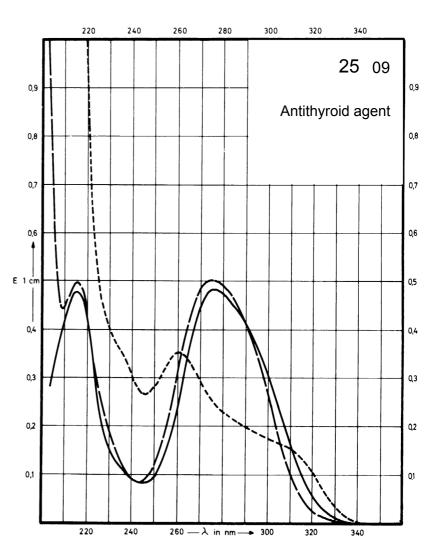

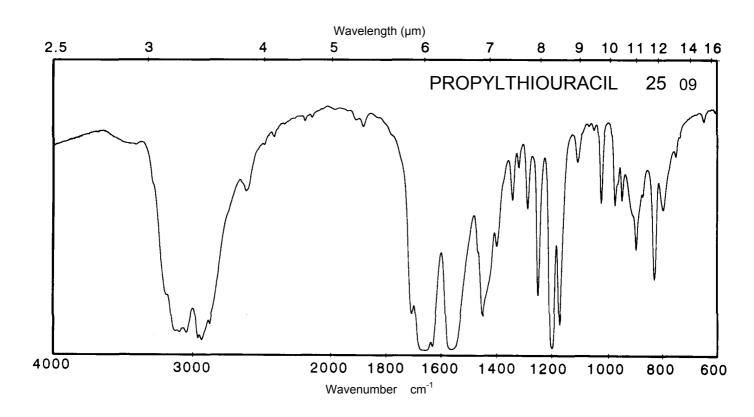





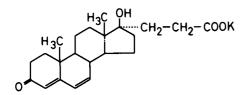


 M_{r} 396.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm		292 nm	294 nm
E ^{1%} 1cm	670		660	654
ε	26570		26180	25940

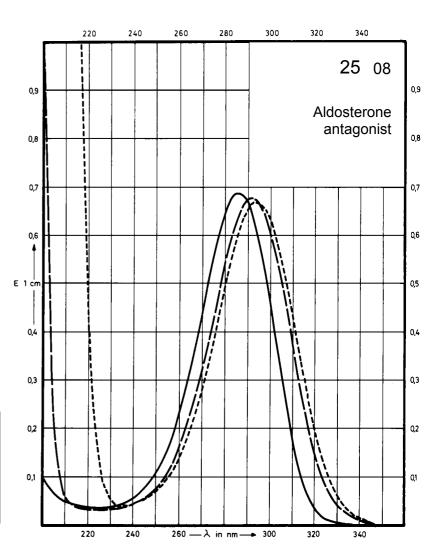

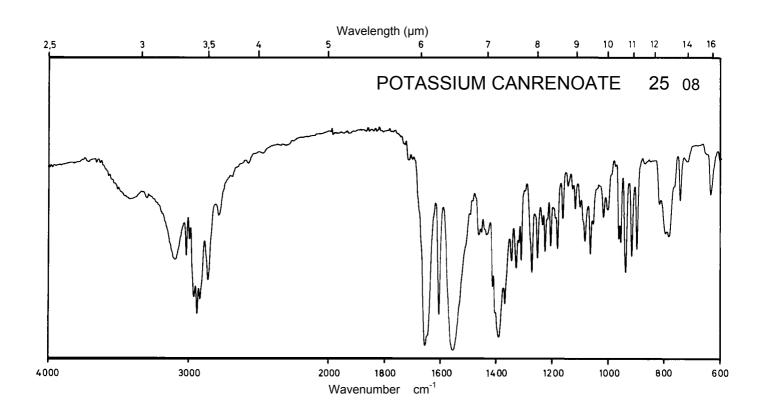




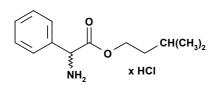


 M_r 257.8


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm			
E ^{1%} _{1cm}	9.2			
3	238			

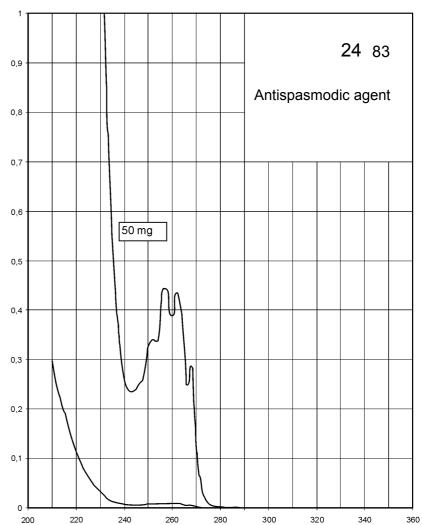

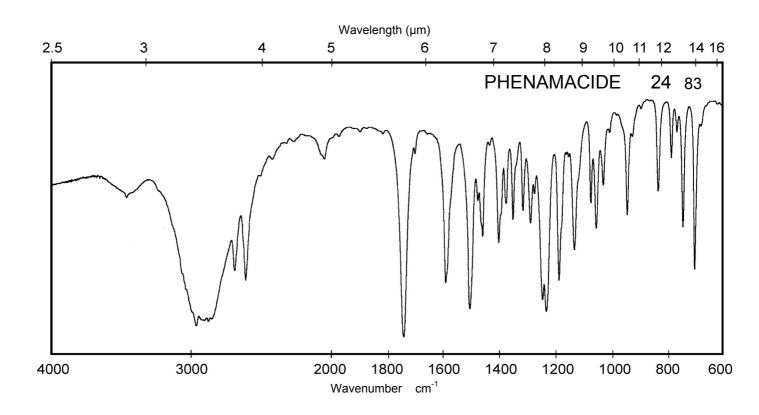





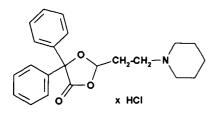



 M_{r} 387.9

Concentration 1 mg / 100 ml 50 mg / 100 ml

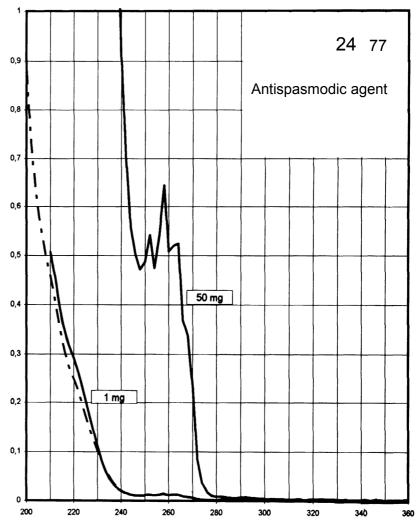


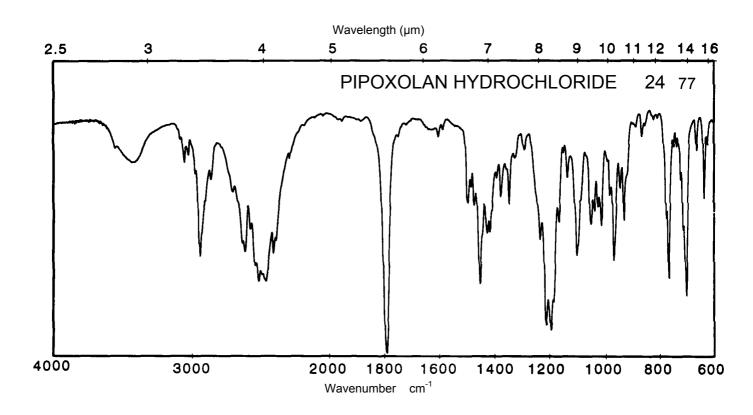

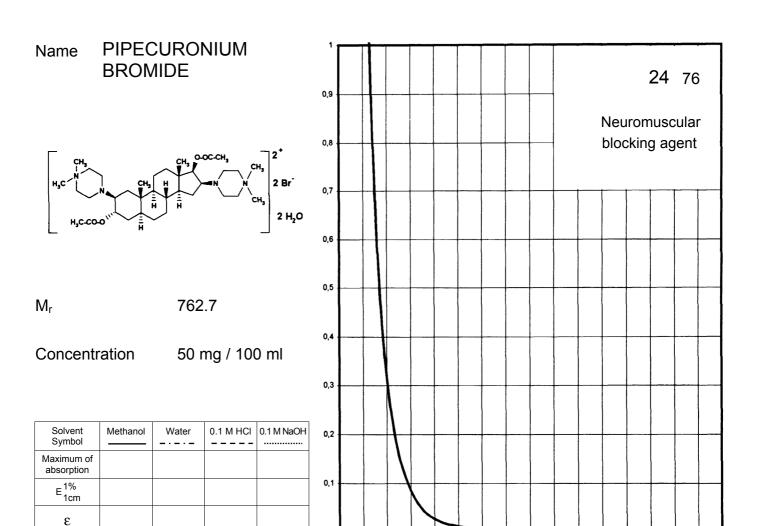

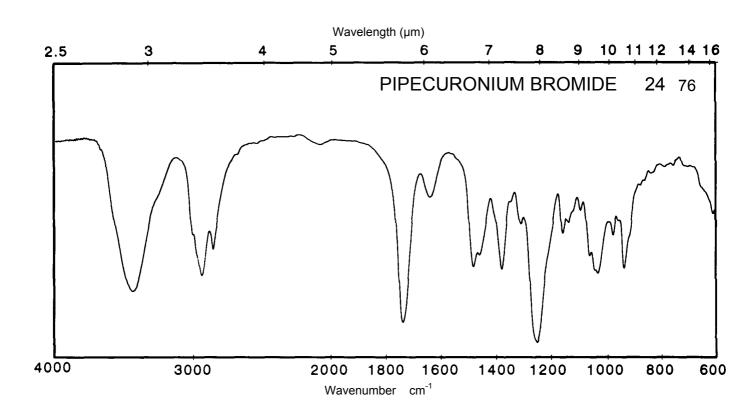





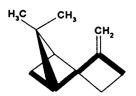


136.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

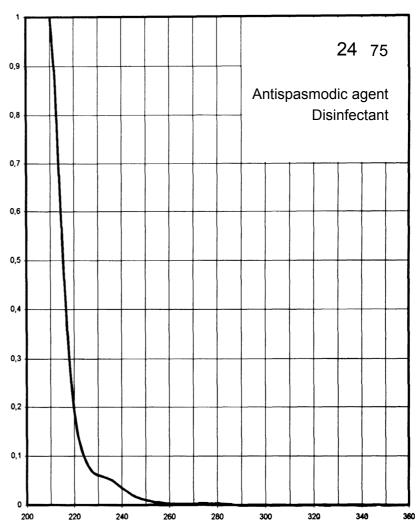

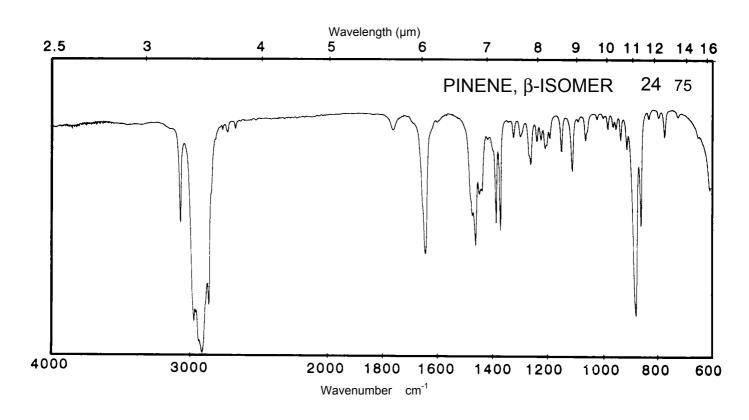




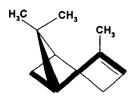


136.2


Concentration 135 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

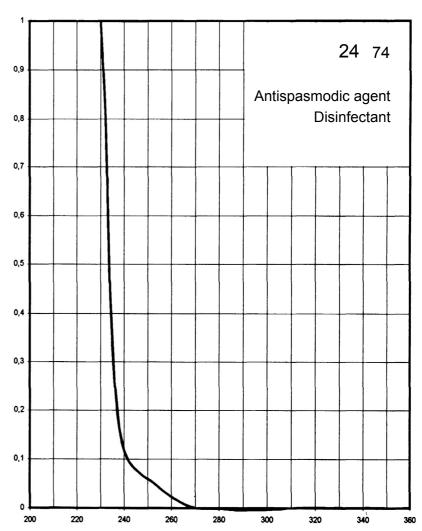

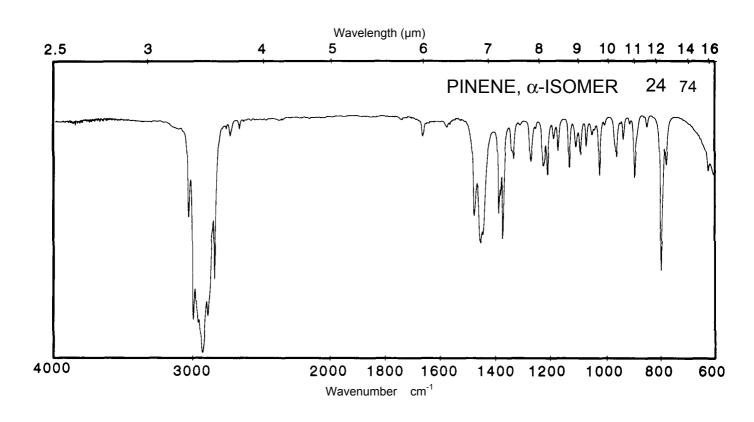




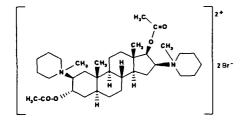


732.7


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

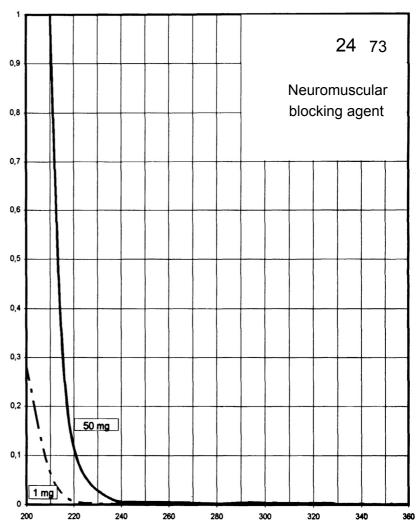

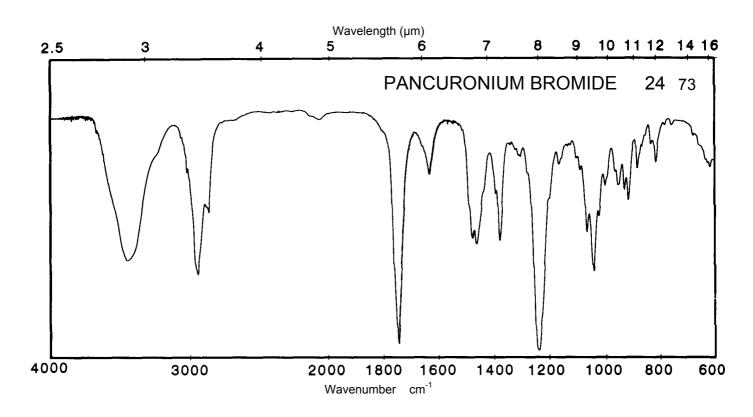




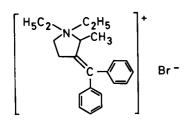


 M_r 386.4


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	244 nm 227 nm		228 nm	
E ^{1%} 1cm	327 347		334	
ε	12600 13400		12900	

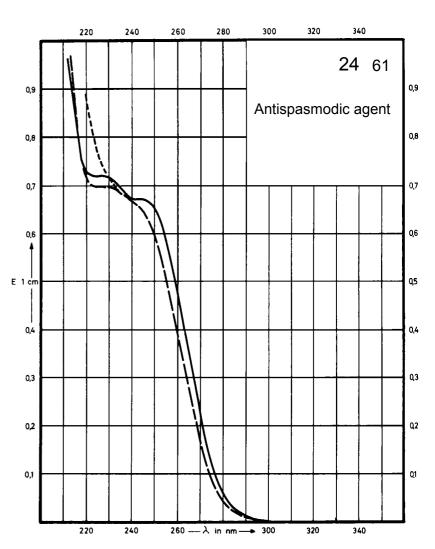

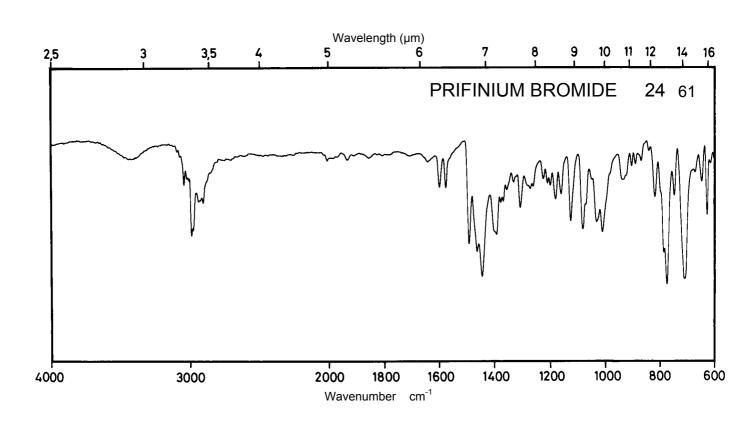




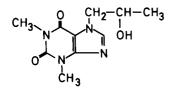


238.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm	273 nm	273 nm	273 nm
E ^{1%} 1cm	376	398	393	400
3	8950	9490	9360	9520

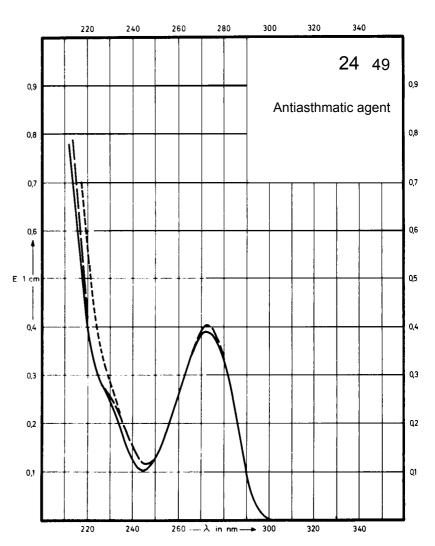

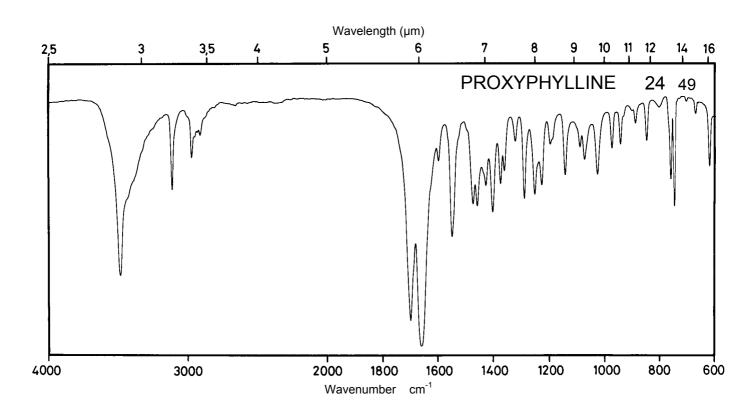

Name PRAMIVERINE HYDROCHLORIDE

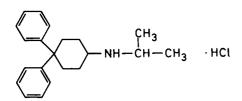



 M_{r}

329.9


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm		261 nm	
E ^{1%} 1cm	14.6		16.3	
ε	480		540	

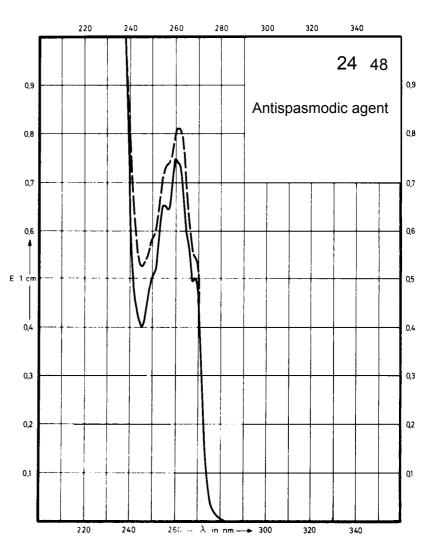

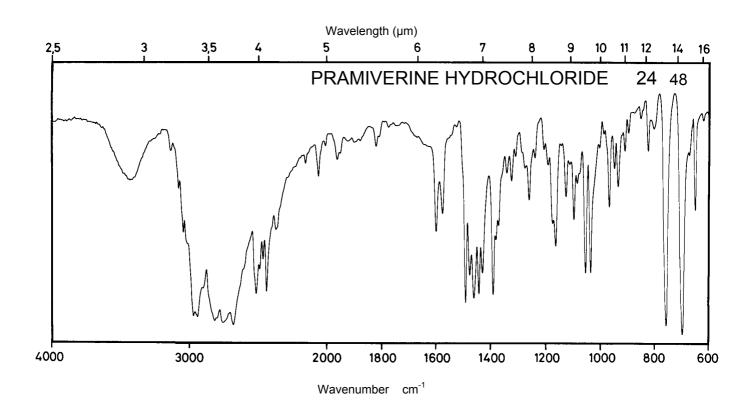




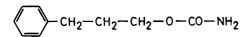


179.2


Concentration 65 mg / 100 ml

	i	i	i	
Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	269 nm			
	259 nm			
absorption	254 nm			
1%	8.8			
E ^{1%} 1cm	11.2			
ICM	9.2			
	157			
3	200			
	165	l	l	l

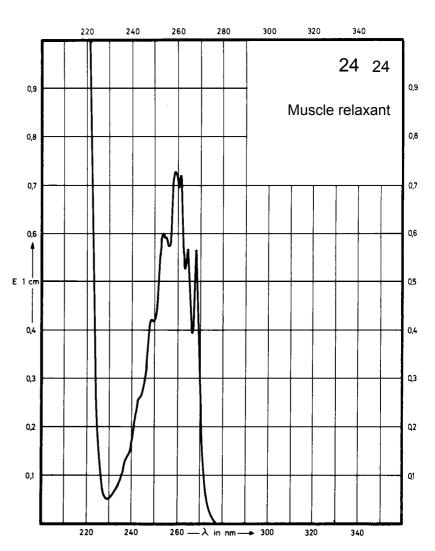

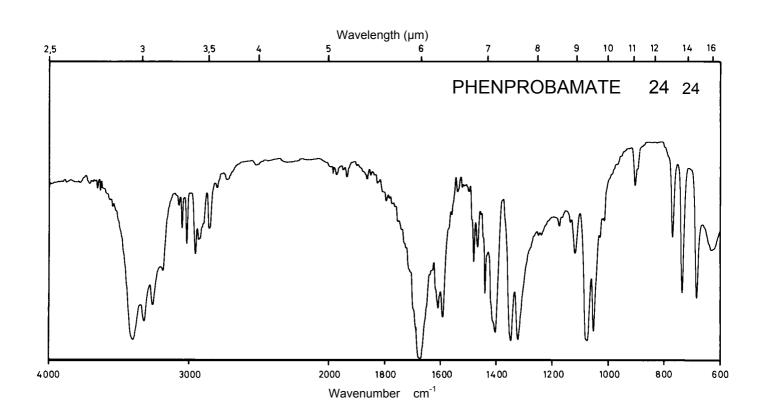

Name PROCYCLIDINE HYDROCHLORIDE

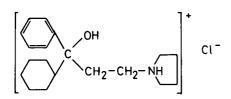



 M_{r}

323.9


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 257 nm 251 nm		263 nm 257 nm 251 nm	
E ^{1%} 1cm	4.50 6.03 4.73		4.50 6.08 4.91	
ε	146 195 153		145 197 159	

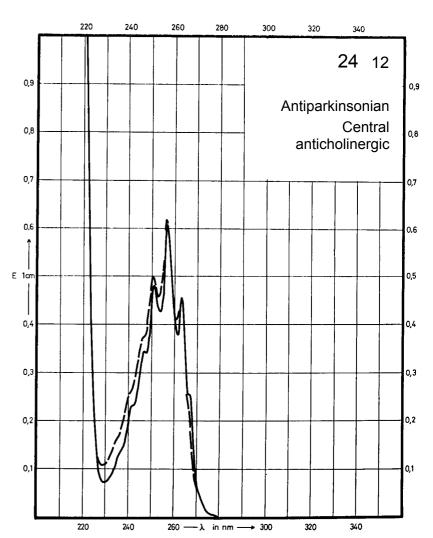

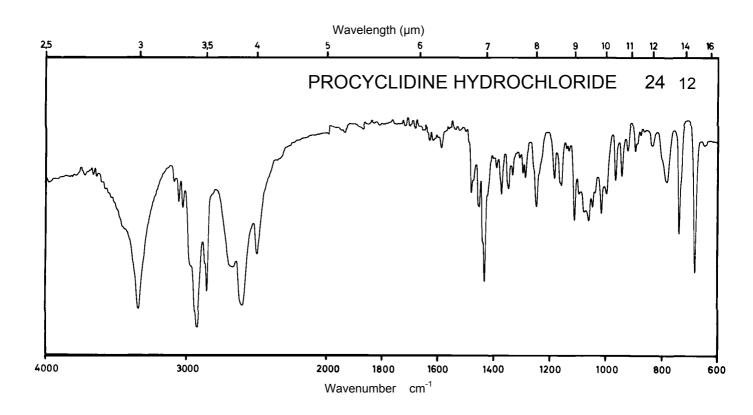




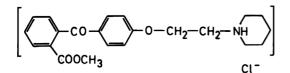


403.9

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm		288 nm	290 nm
E ^{1%} 1cm	420		450	430
ε	16960		18180	17370

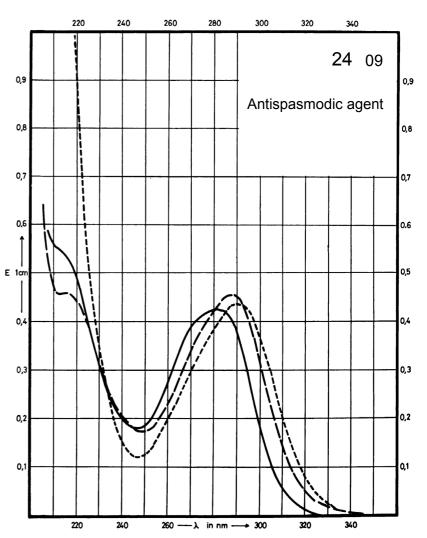

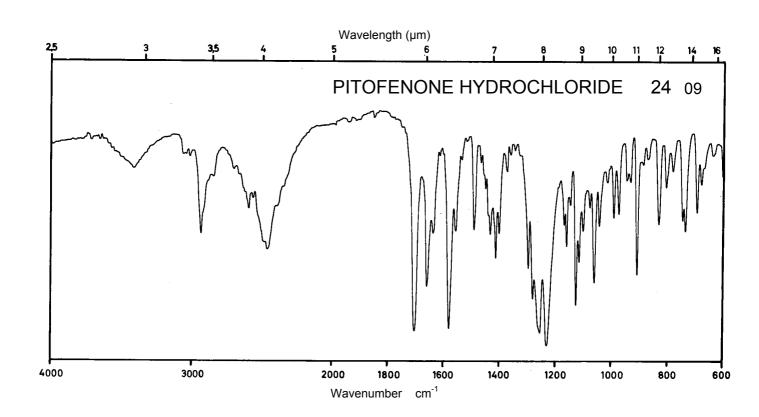

Name PAPAVERINE HYDROCHLORIDE


 M_{r}


375.9

Concentration 0.5 mg / 100 ml 2.0 mg / 100 ml


Solvent Symbol	Methanol*	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm 250 nm		308 nm 284 nm 249 nm	326 nm 277 nm 236 nm
E ^{1%} 1cm	244 1910		228 178 1700	117 196 1800
ε	9170 71800		8570 6690 63900	4400 7370 67670

* 1 M HCl + Methanol (1 + 99)

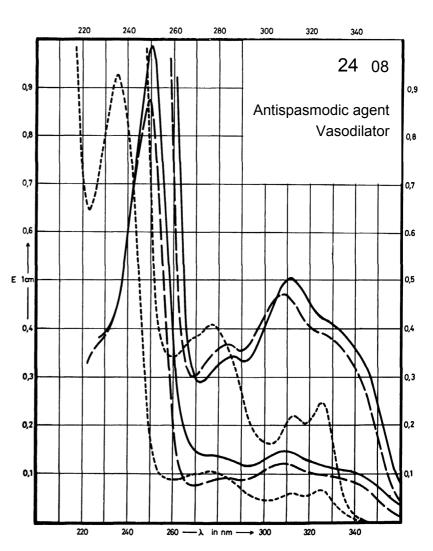

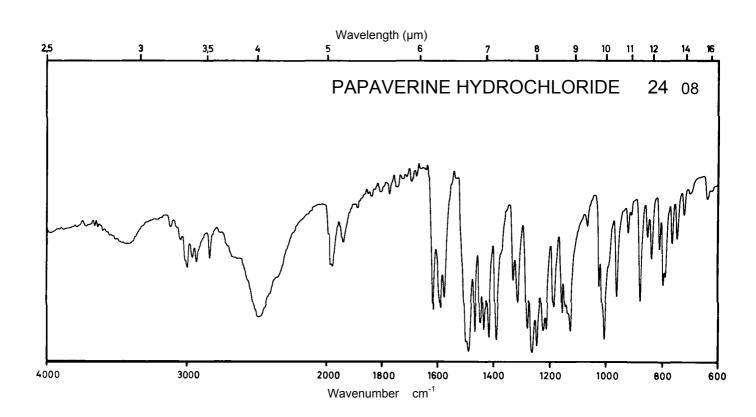



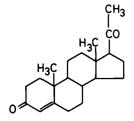



314.5

Concentration


0.9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm			
E ^{1%} 1cm	540			
ε	16980			

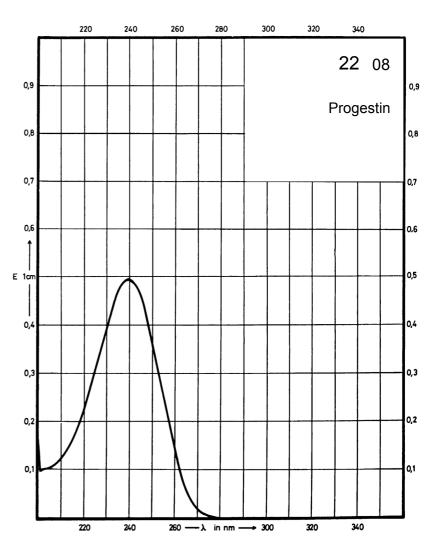

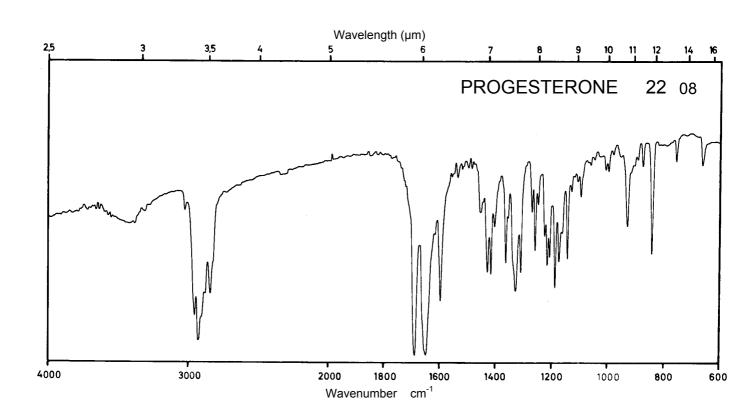

Name PAROXETINE HYDROCHLORIDE

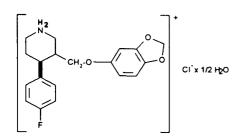



 M_{r}

374.8


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	294 nm 235 nm	293 nm 233 nm	293 nm 233 nm	293 nm
E ^{1%} 1cm	125 121	110 112	110 112	101
3	4700 4500	4100 4200	4100 4200	3800

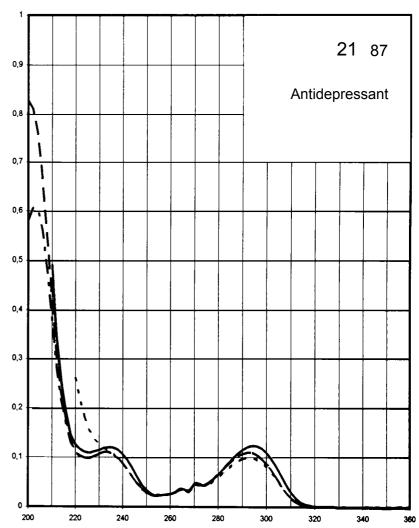

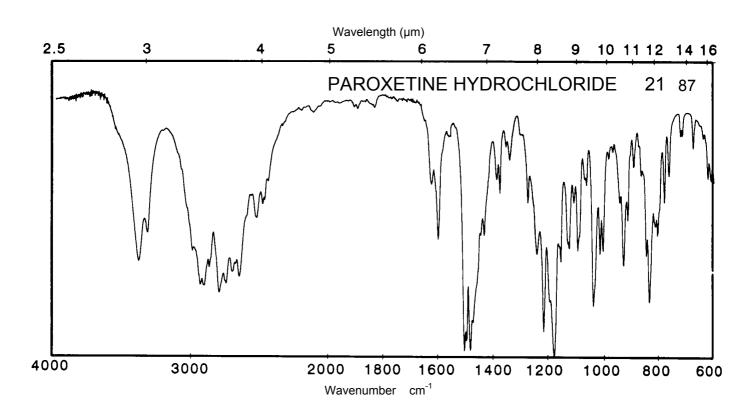




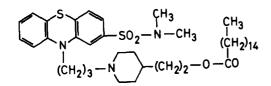


714.1

Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	313 nm 265 nm 236 nm			
E ^{1%} 1cm	51 498 265			
ε	3600 35600			

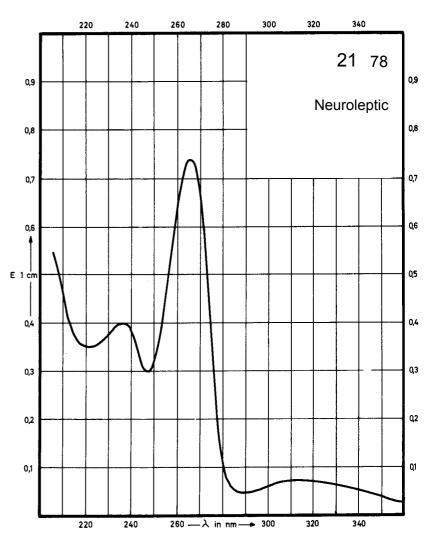

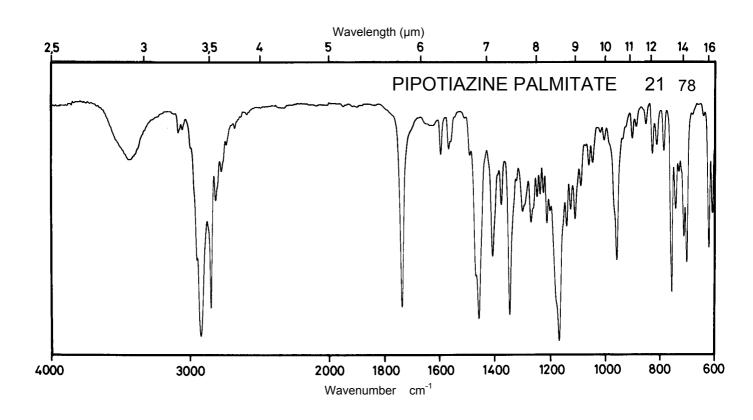



Name PIRISUDANOL DIMALEATE

 M_r 572.5


Concentration 4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm	326 nm 291 nm	291 nm	309 nm
E ^{1%} 1cm	117	64 94	159	121
ε	6700	3660 5380	9110	6950

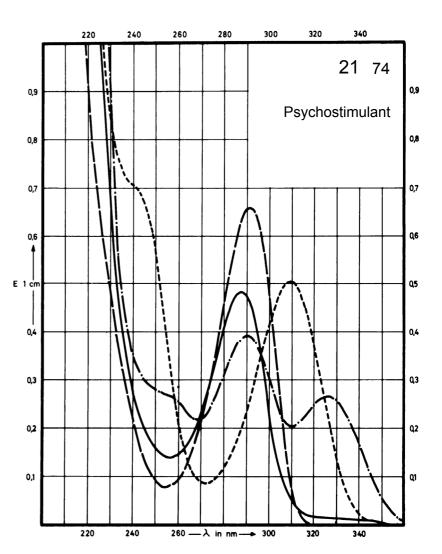

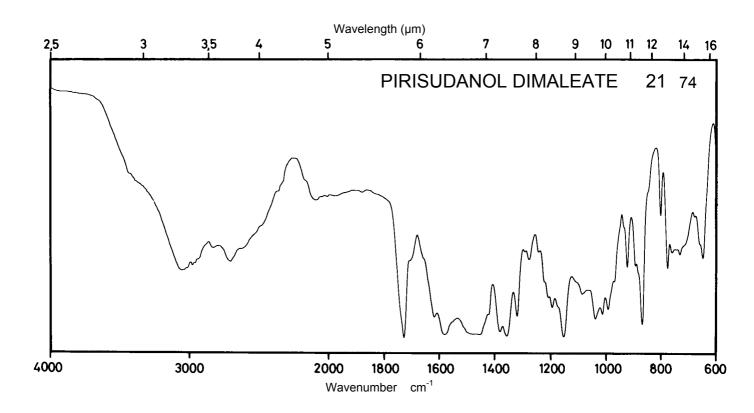


Name PIPRADOL HYDROCHLORIDE

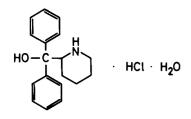


 M_r 321.9


Concentration 30 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm 254 nm		259 nm	
E ^{1%} 1cm	14.5 12.2		14.8	
3	470 390		475	

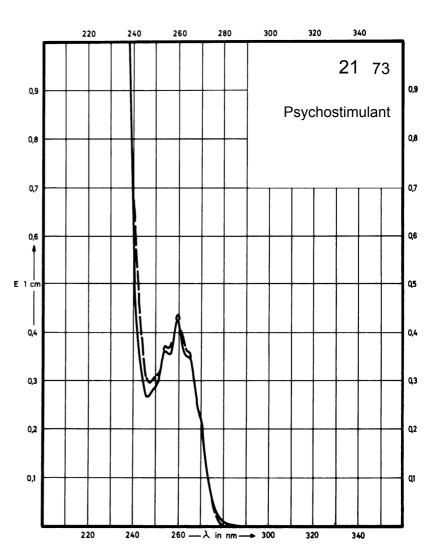

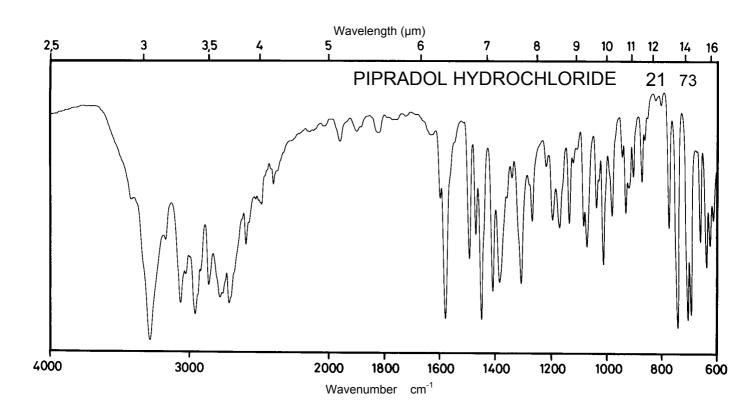




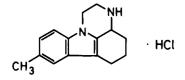


262.8


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm	276 nm	276 nm	289 nm
E 1% 1cm	303	269	271	258
ε	7960	7070	7120	6780

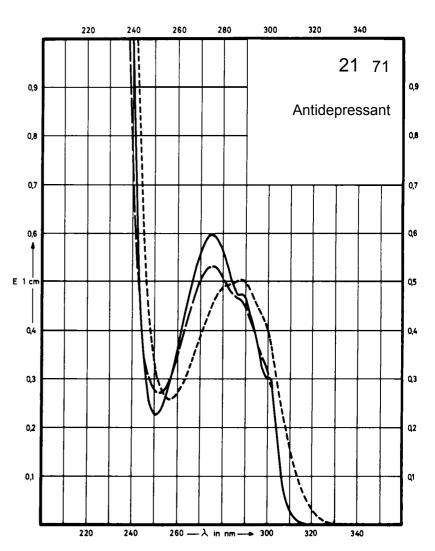

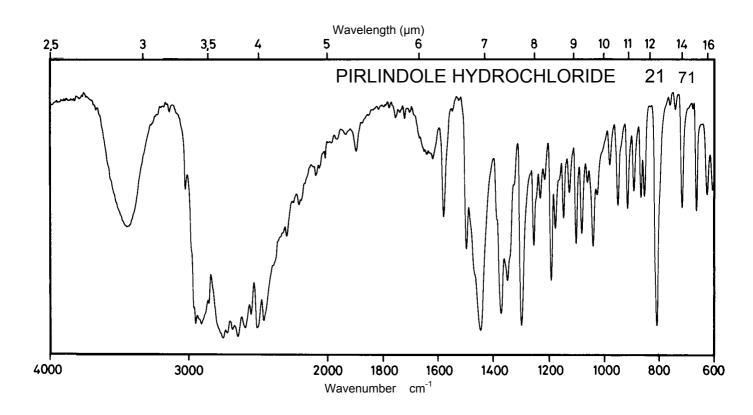




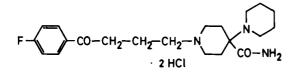


 M_r 448.4


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	244 nm		248 nm	247 nm
E ^{1%} 1cm	287		273	274
3	12900		12200	12300

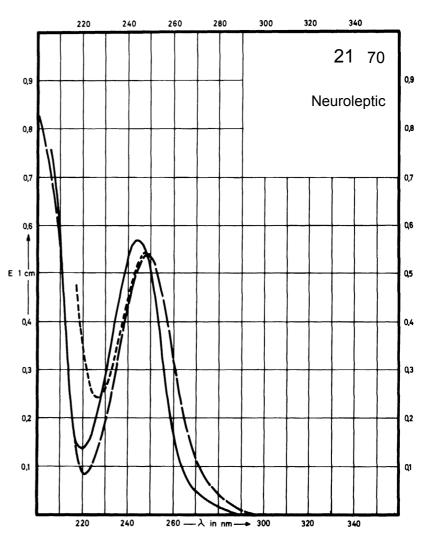

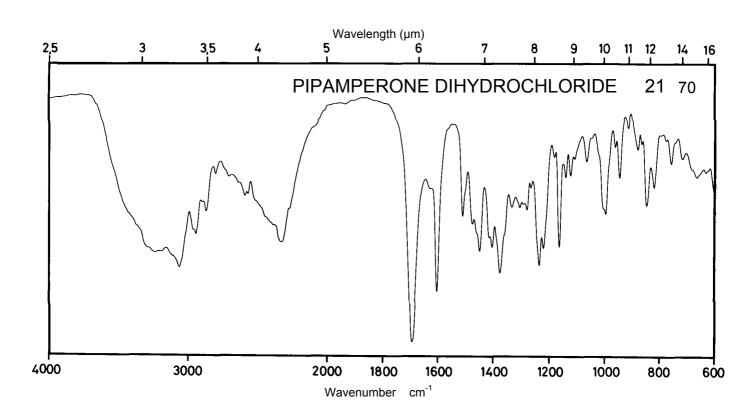



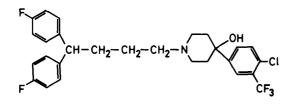



 M_r 524.0

Concentration 10 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI*	0.1 M NaOH
Maximum of absorption	273 nm 267 nm		273 nm 267 nm	
E ^{1%} 1cm	54 52		48 47.5	
ε	2800 2700		2520 2490	

^{*} cont. 40% methanol

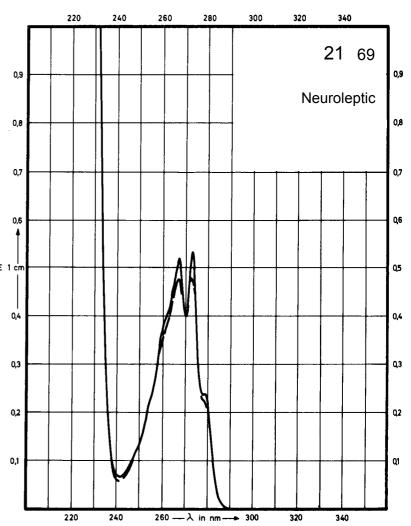

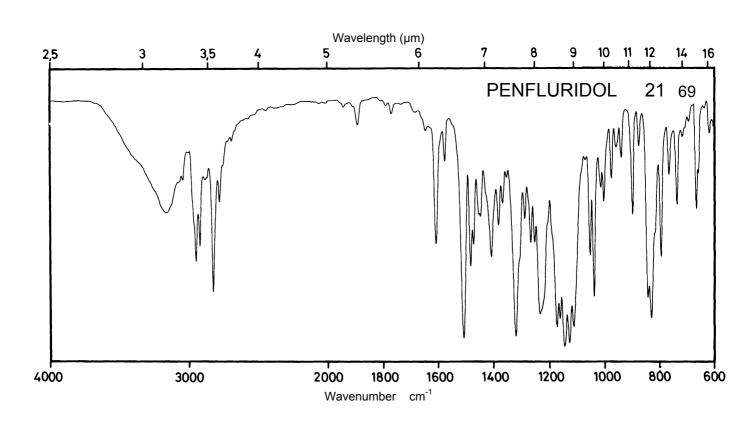




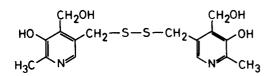


368.5

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	299 nm		297 nm	315 nm
E ^{1%} 1cm	515		490	421
ε	18970		18050	15500

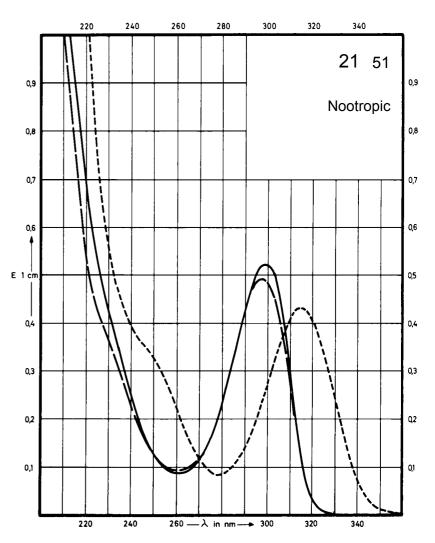

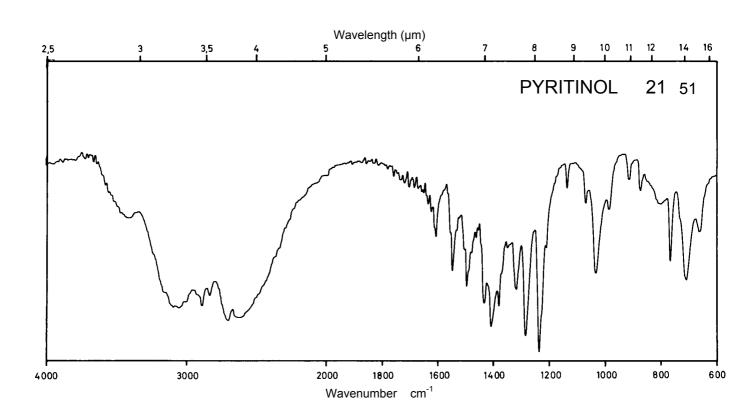



Name PROTRIPTYLINE HYDROCHLORIDE

 M_r 299.8


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	293 nm		292 nm	292 nm
E ^{1%} 1cm	466		451	425
ε	13960		13520	12750

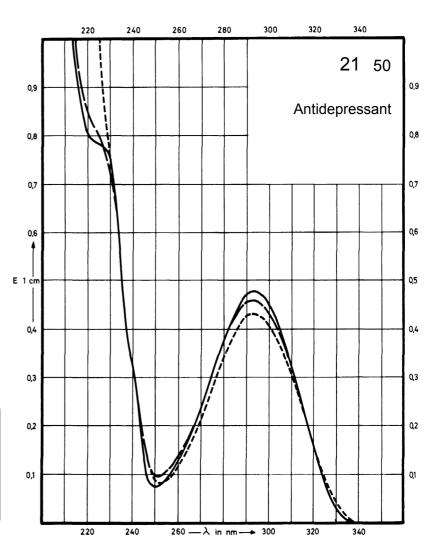

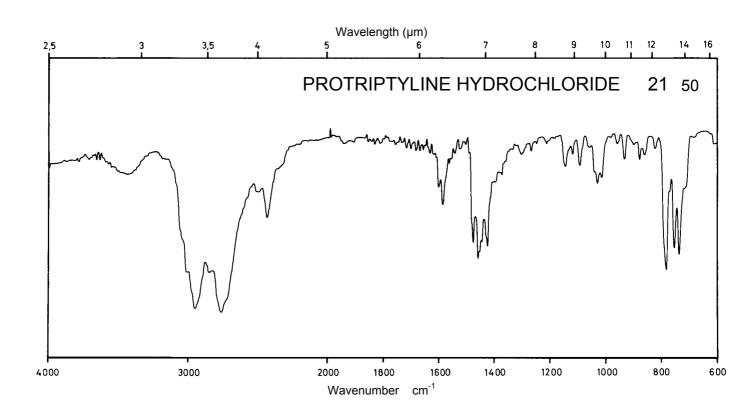

Name PHENDIMETRAZINE TARTRATE

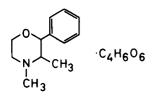



 M_r

341.4


Concentration 73 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 257 nm 251 nm		261 nm 257 nm 251 nm	263 nm 257 nm 251 nm
E ^{1%} 1cm	5.2 6.2 4.6		6.1 6.8 5.0	5.2 6.4 4.9
ε	176 212 157		210 230 170	176 218 168

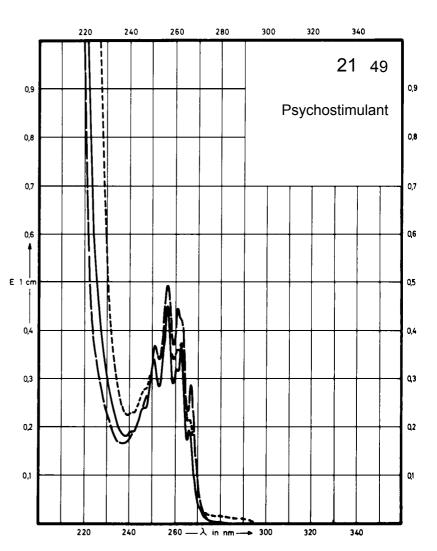

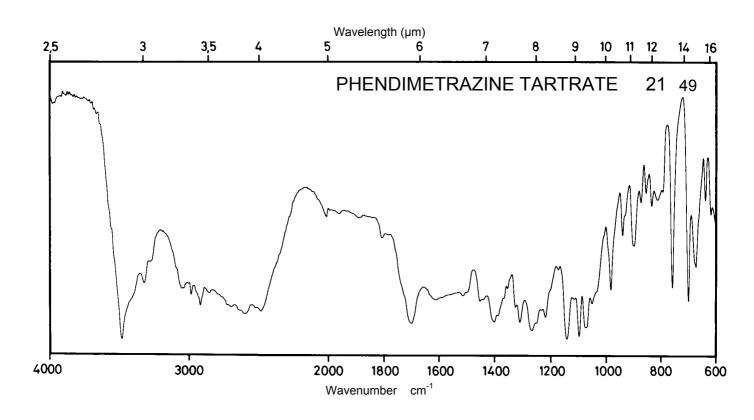




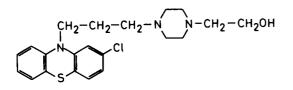


 M_r 404.0

Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	312 nm 258 nm		307 nm 255 nm	311 nm 257 nm
E ^{1%} 1cm	110 895		100 800	101 783
ε	4450 36140		4020 32350	4080 31640

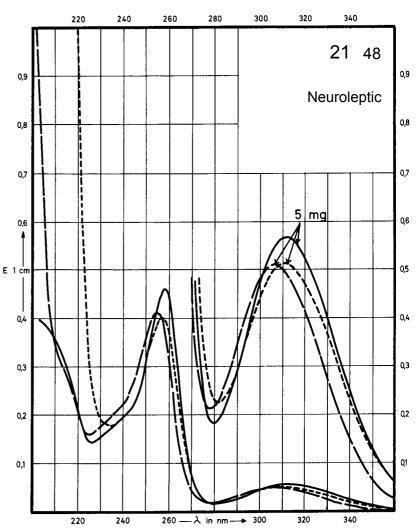

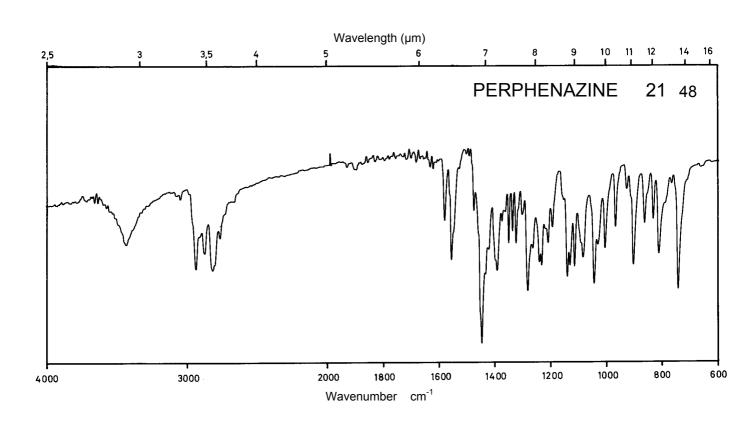

Name PROMAZINE PHOSPHATE


$$\begin{array}{c|c} \text{S} & \cdot \text{H}_3\text{PO}_4 \\ \vdots \\ \text{CH}_2\text{-CH}_2\text{-CH}_2\text{-N} \\ \text{CH}_3 \end{array}$$


 M_{r}

382.4

Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	303 nm 253 nm		301 nm 252 nm	304 nm 253 nm
E ^{1%} 1cm	105 860		101 800	97 750
ε	4020 32890		3860 30590	3710 28680

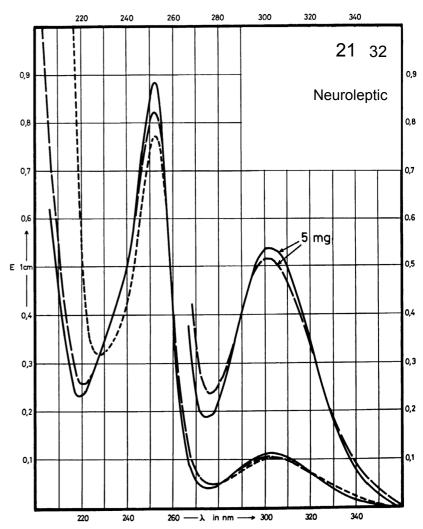

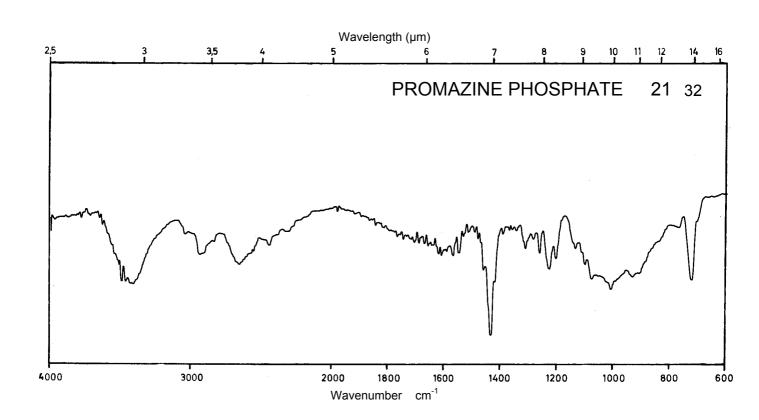


Name PROCHLORPERAZINE MALEATE


 M_{r}

606.1

Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm 257 nm		307 nm 254 nm	310 nm 256 nm
E ^{1%} 1cm	73 620		67 560	65 500
ε	4420 37580		4060 33940	3940 30300

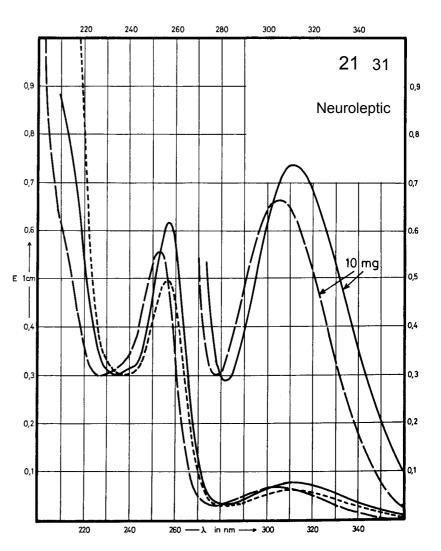

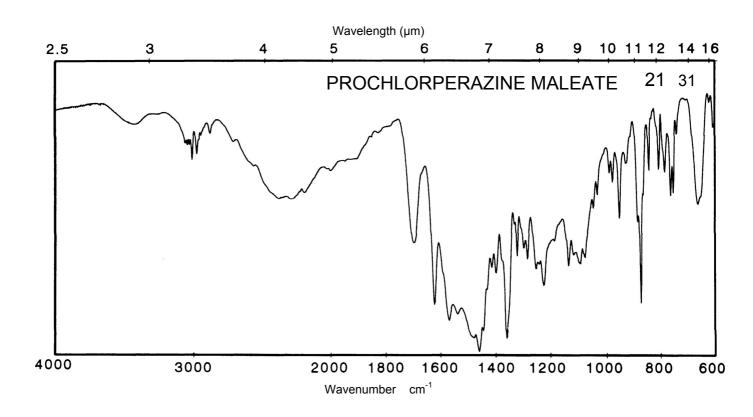


Name PIMOZIDE


 M_{r}

461.6


Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm		278 nm	
E ^{1%} 1cm	145		145	
3	6690		6690	

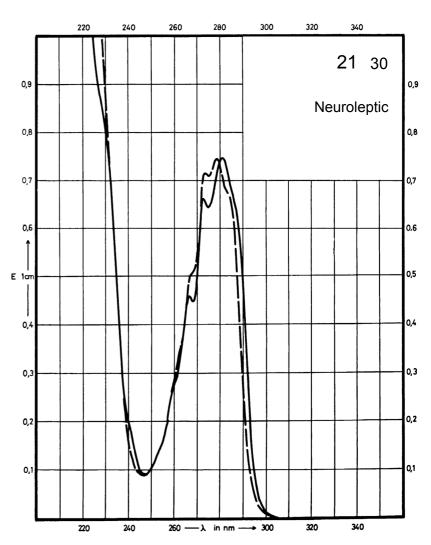

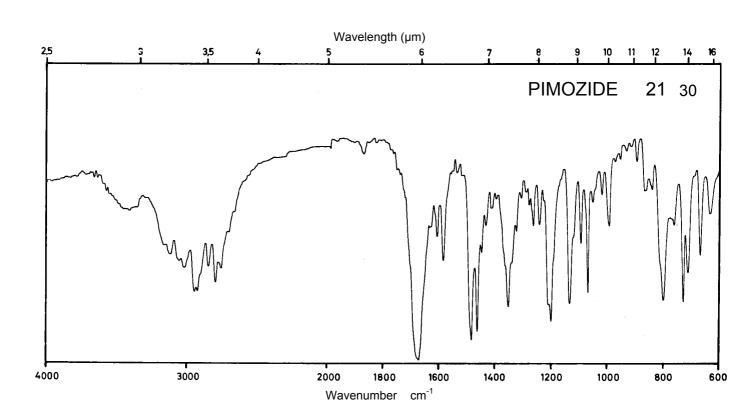




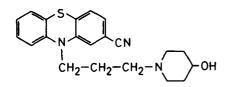


365.5

Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm 271 nm 232 nm		268 nm 232 nm	270 nm
E ^{1%} 1cm	77 900 596		763 648	767
ε	2810 32900 21780		27890 23680	28030

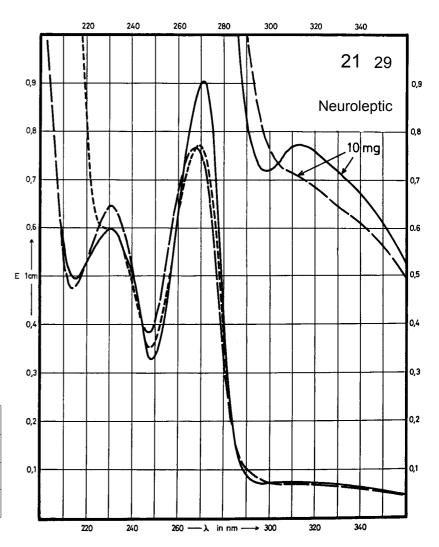

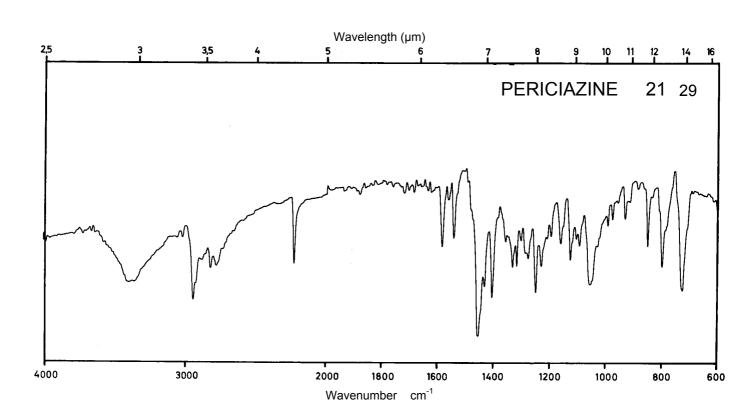


Name PERAZINE DIMALONATE


$$\begin{array}{c|c}
S & & \\
N & & \\
CH_2-CH_2-CH_2-N & N-CH_3 \\
& \cdot 2 & C_3H_4O_4
\end{array}$$

 M_r 547.6


Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	307 nm 255 nm		302 nm 252 nm	307 nm 254 nm
E ^{1%} 1cm	79 616		71 560	75 564
ε	4520 35210		4060 32010	4290 32240

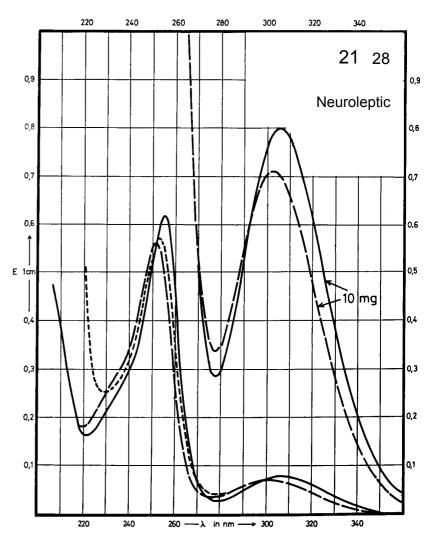

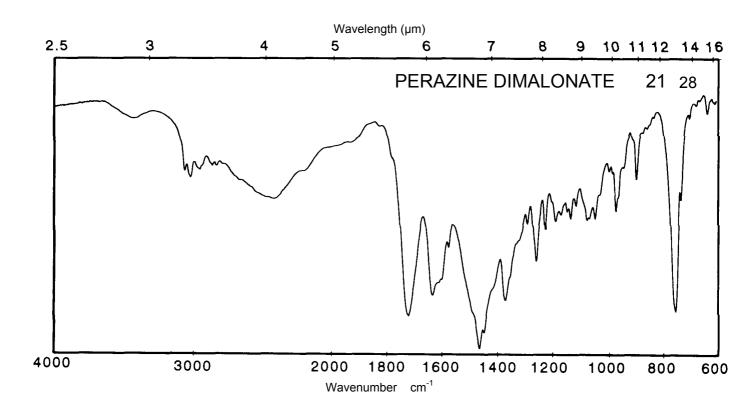



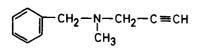


· HCI


 M_{r}

195.7


Concentration 60 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm 262 nm 256 nm		268 nm 262 nm 256 nm	263 nm 257 nm 252 nm
E ^{1%} 1cm	11.1 15.1 13.9		11.2 15.5 14.2	7.8 10.2 8.3
ε	220 300 270		220 300 280	150 200 160

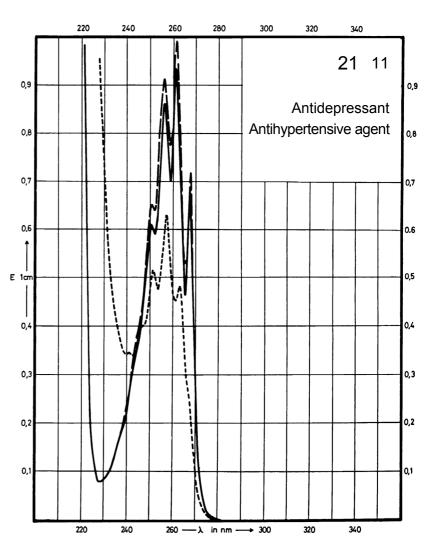

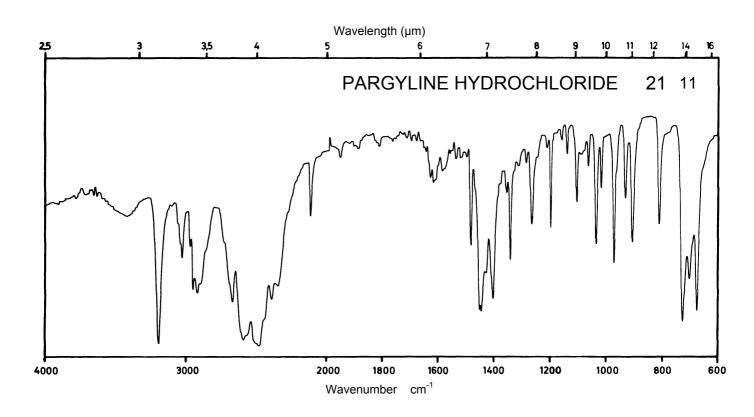




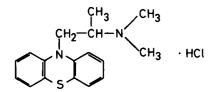


320.9


Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	300 nm 252 nm		298 nm 248 nm	254 nm
E 1% 1cm	109 922		109 929	949
ε	3500 29590		3500 29810	30450

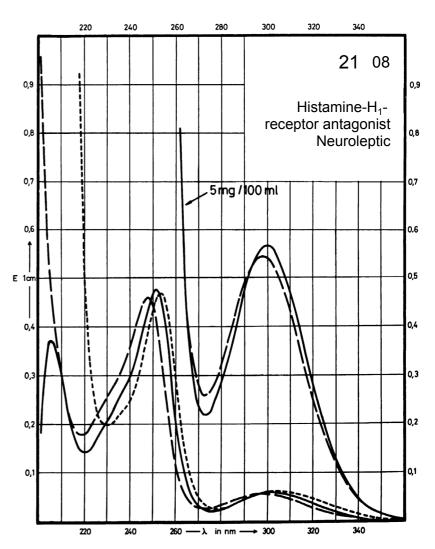

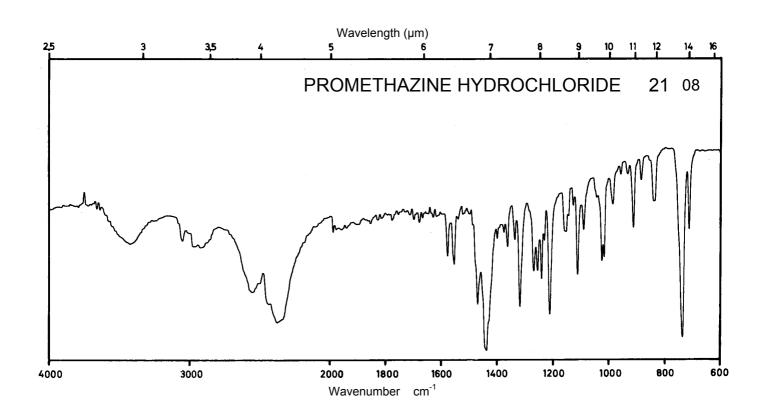

Name PROFENAMINE HYDROCHLORIDE

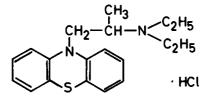


 M_{r}


348.9

Concentration

1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	300 nm 250 nm		300 nm 250 nm	
E ^{1%} 1cm	92 820		88 798	
ε	3210 28610		3070 27840	

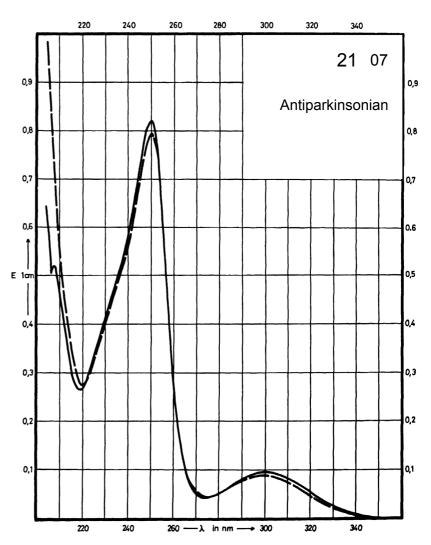

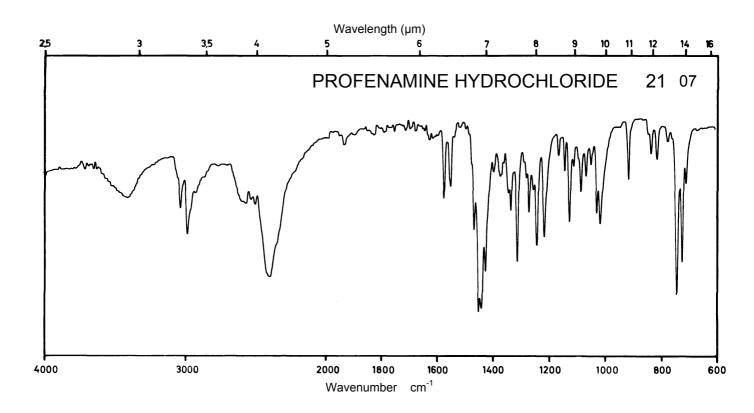

Name PRENALTEROL HYDROCHLORIDE


$$\begin{array}{c} \text{OH} \\ \text{I} \\ \text{O-CH}_2\text{-CH-CH}_2\text{-NH-CH} \\ \\ \cdot \text{HCI} \end{array}$$


 M_{r}

261.8


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm	287 nm	287 nm	305 nm
	225 nm	222 nm	222 nm	235 nm
E ^{1%}	104	93	93	109
1cm	300	285	285	382
ε	2720	2430	2430	2850
	7860	7460	7460	10000

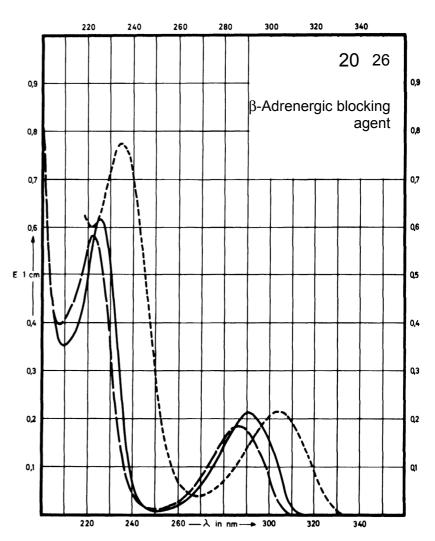

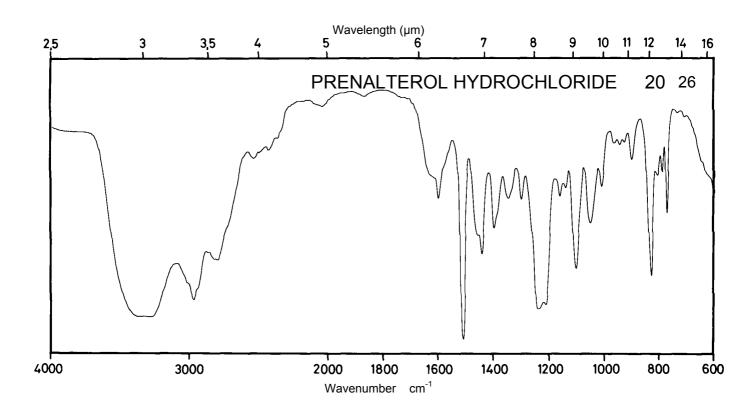

Name PINDOLOL

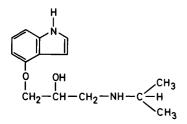


 M_{r}


248.3

Concentration


1.7 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 265 nm		287 nm 264 nm	285 nm 262 nm
E ^{1%} 1cm	169 320		153 292	158 299
ε	4200 7950		3800 7250	3920 7420

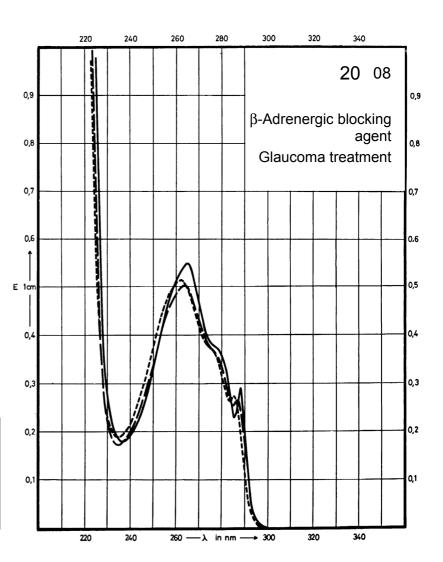

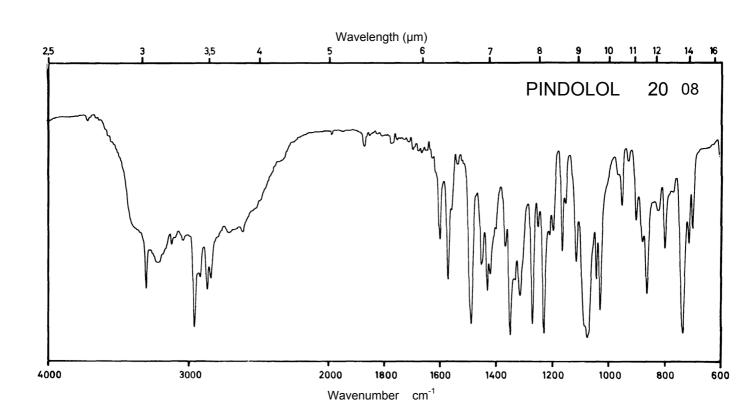




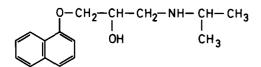


259.4

Concentration 2.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	320 nm 292 nm			
E ^{1%} 1cm	80 236			
ε	2070 6110			

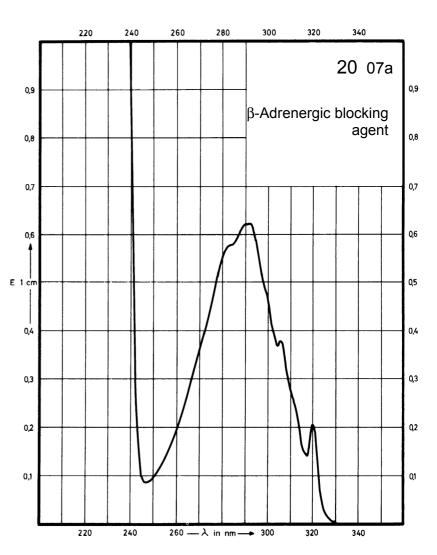

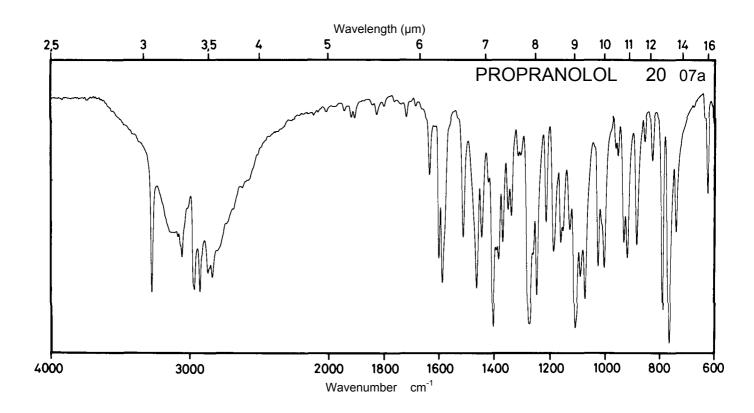



Name PROPRANOLOL HYDROCHLORIDE

 M_r 295.8

Concentration 3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	319 nm 289 nm		319 nm 288 nm	319 nm 291 nm
E ^{1%} 1cm	66 207		68 196	75 199
ε	1950 6120		2010 5800	2220 5890

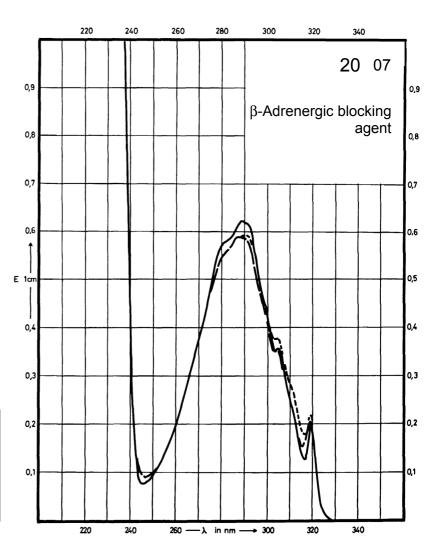

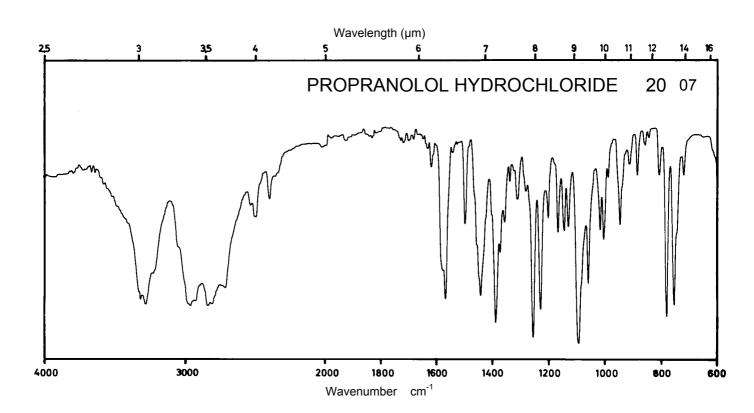


Name PRACTOLOL


 M_{r}

266.3

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	248 nm		243 nm	244 nm
E ^{1%} 1cm	610		450	456
ε	16240		12000	12140

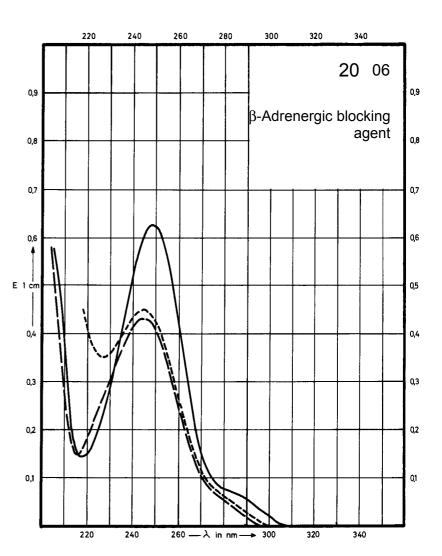

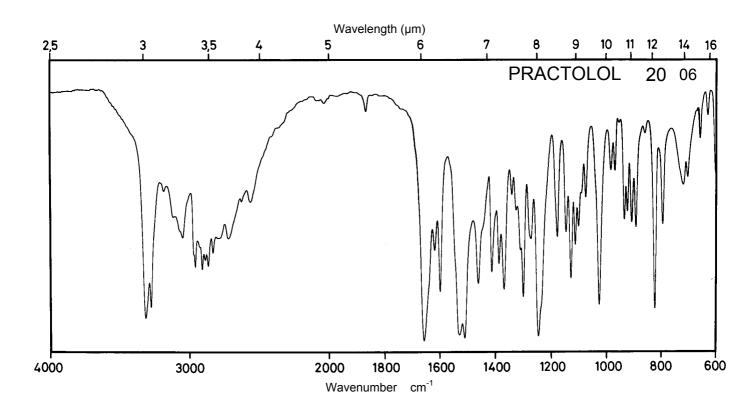


Name PYRIDOSTIGMINE BROMIDE


 M_r

261.1

Concentration 2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm		270 nm	Decom- position observed
E ^{1%} _{1cm}	185		182	
ε	4830		4760	

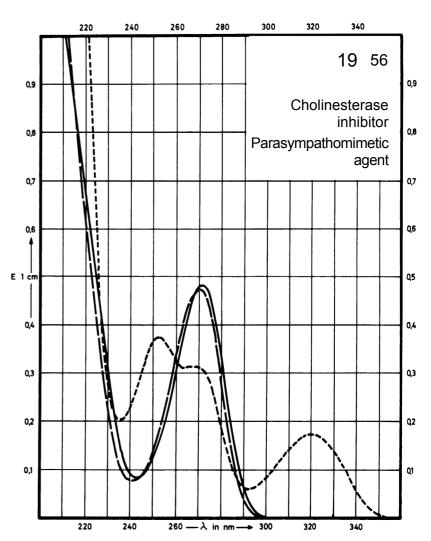

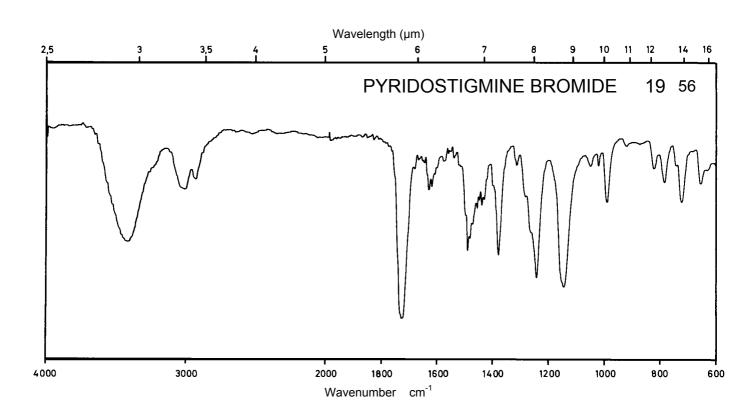

Name PHYSOSTIGMINE


 M_{r}


275.3

Concentration


1.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm 252 nm		302 nm 245 nm	298 nm
E ^{1%} 1cm	109 473		103 463	537
3	3000 13020		2840 12750	14780

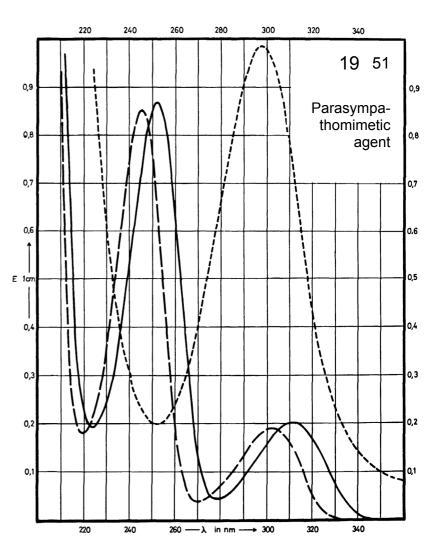

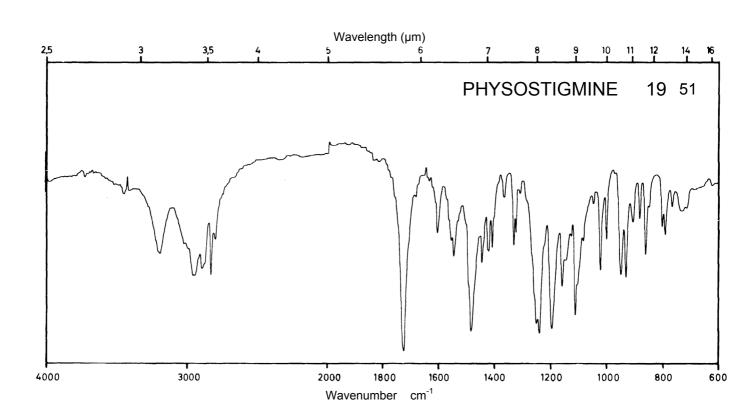




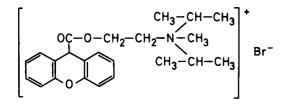


 $\,M_r\,$


448.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm	281 nm	281 nm	282 nm
	246 nm	242 nm	242 nm	249 nm
E ^{1%}	61	60	60	58
1cm	118	105	105	123
ε	2760	2680	2680	2610
	5270	4690	4690	5520

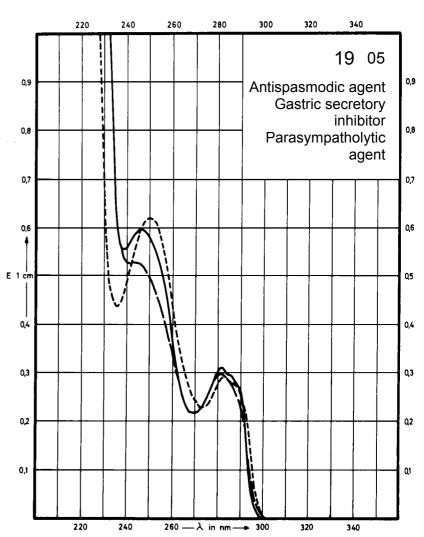

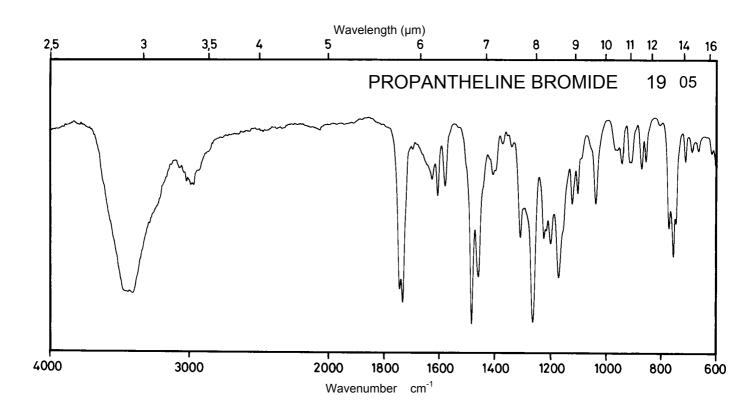





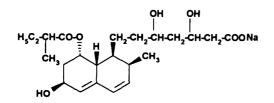


446.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	237 nm	239 nm	238 nm	238 nm
E ^{1%} 1cm	510	469	451	473
3	22700	20900	20100	21100

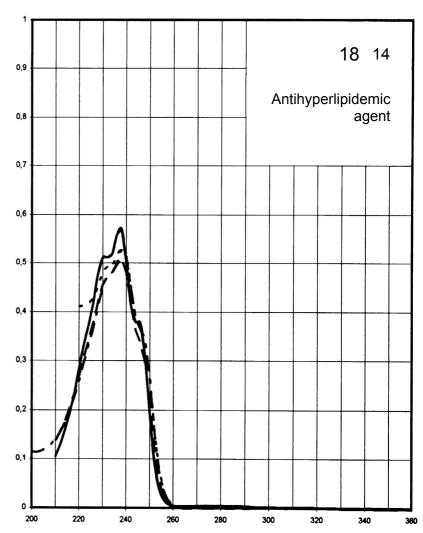

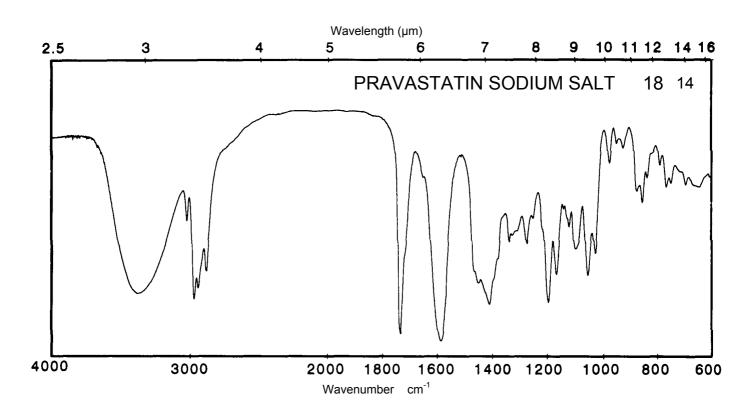





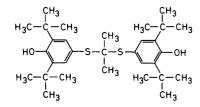



 M_r 516.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm			
E ^{1%} 1cm	428			
ε	22140			

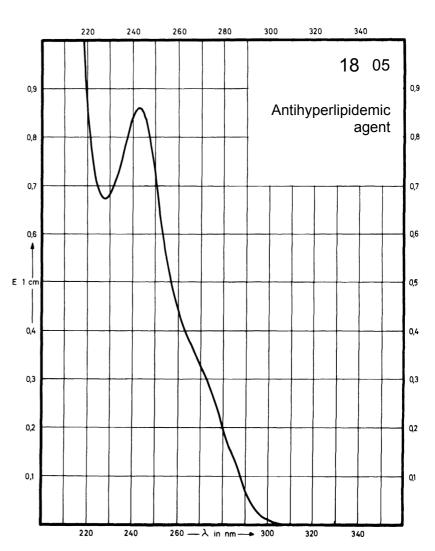

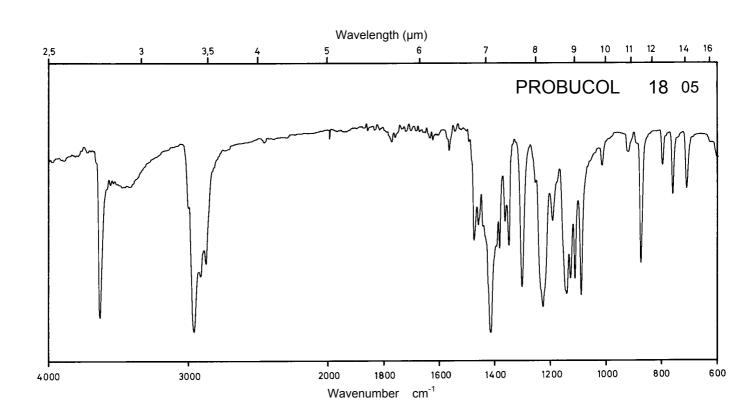

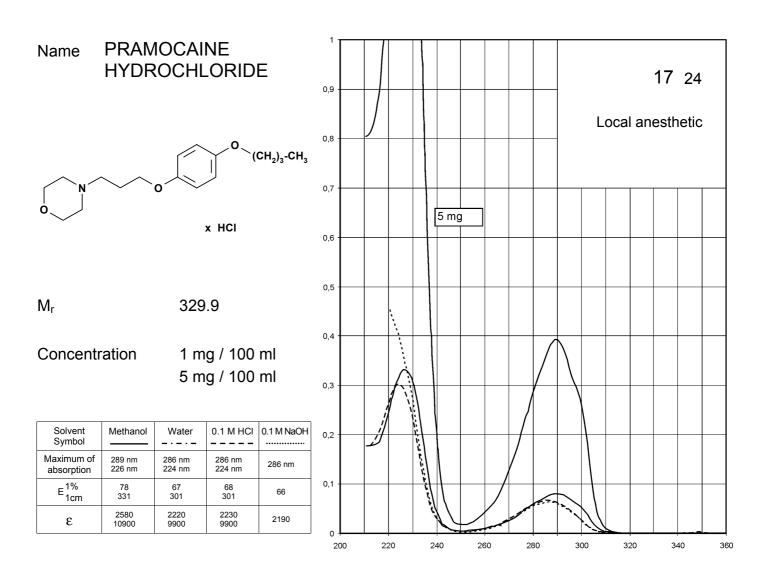

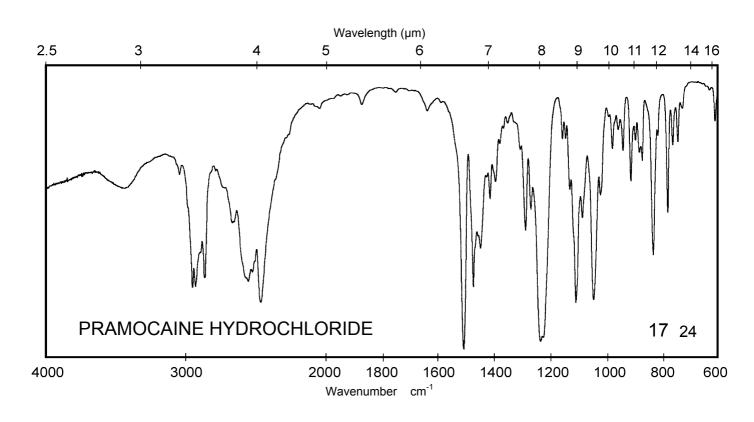



Name PROXYMETACAINE HYDROCHLORIDE

 M_r


330.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	313 nm 266 nm		257 nm	Decom- position observed
E ^{1%} 1cm	145 276		560	
ε	4800 9100		18500	

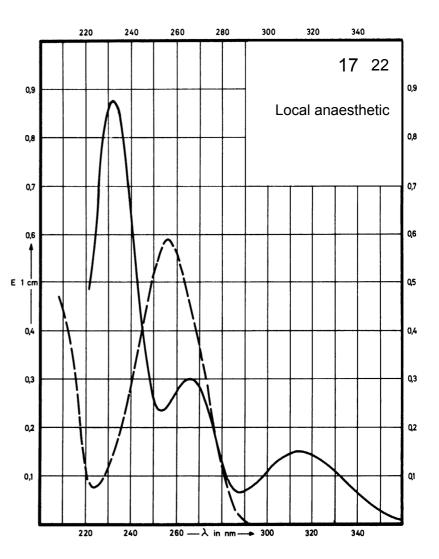

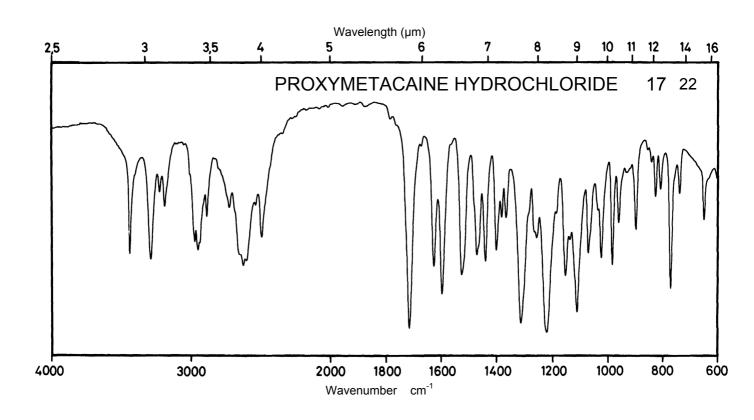





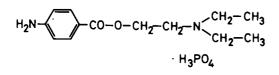


 M_r 334.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	297 nm 223 nm	290 nm 220 nm	280 nm 227 nm	
E ^{1%} 1cm	672 252	549 256	53 385	
ε	22500 8440	18400 8550	1770 12900	

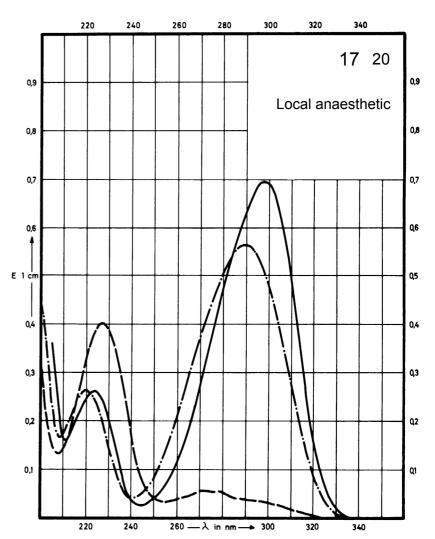

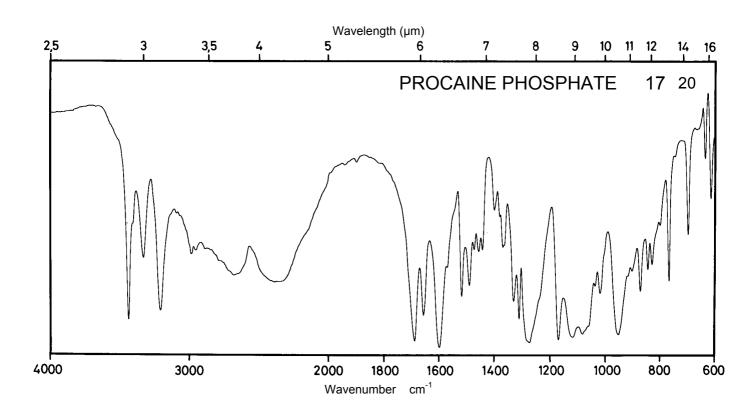





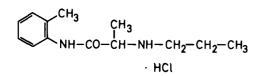


 M_r 256.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	229 nm			
E ^{1%} 1cm	260			
3	6670			

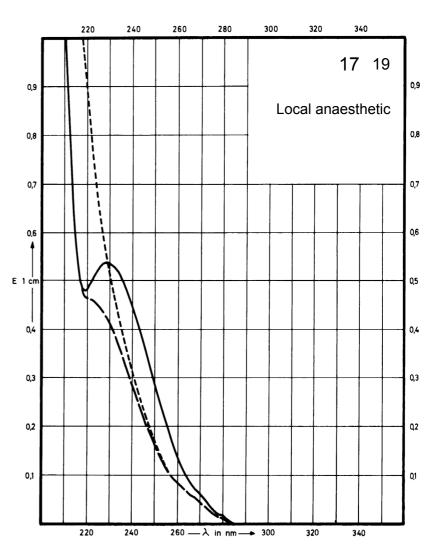

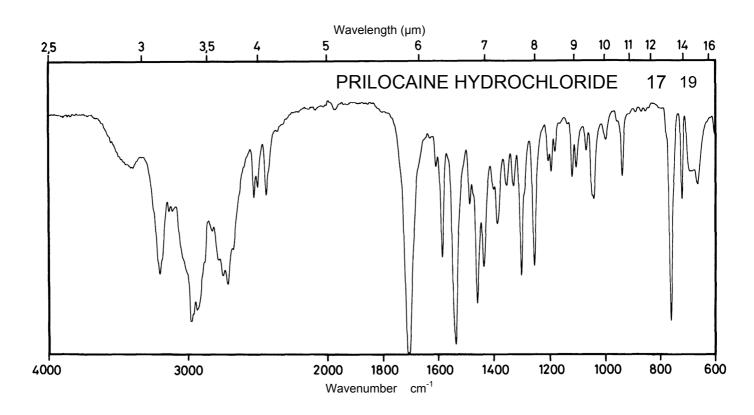


Name POLYCAINE HYDROCHLORIDE

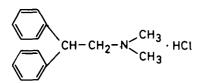


 M_{r}

261.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm	258 nm	258 nm	
E ^{1%}	15.9	16.7	16.9	
ε	416	437	442	

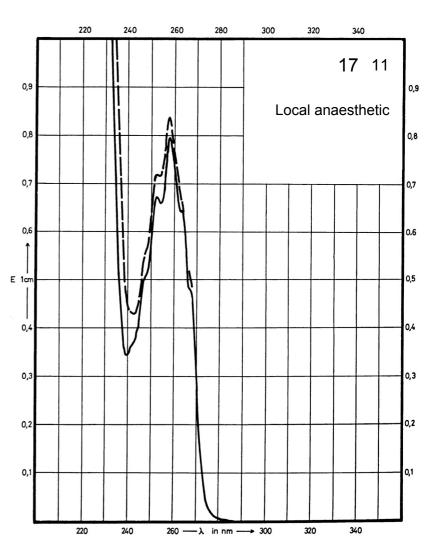

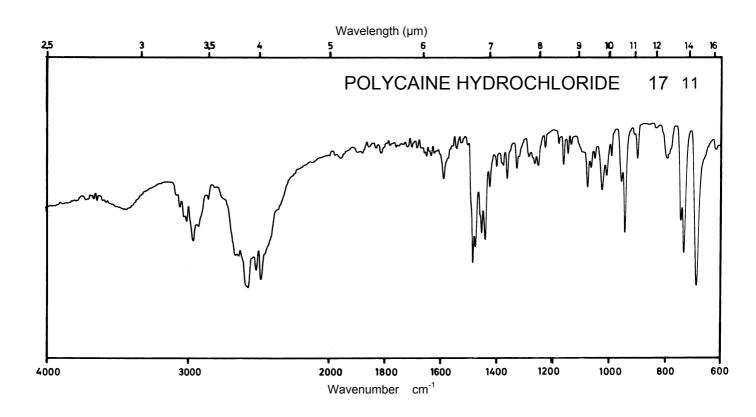



Name PROCAINE HYDROCHLORIDE

 M_{r}


272.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	296 nm	290 nm 220 nm	280 nm 227 nm	268 nm
E ^{1%} 1cm	815	684 328	94 529	447
ε	22230	18660 8950	2560 14430	12190

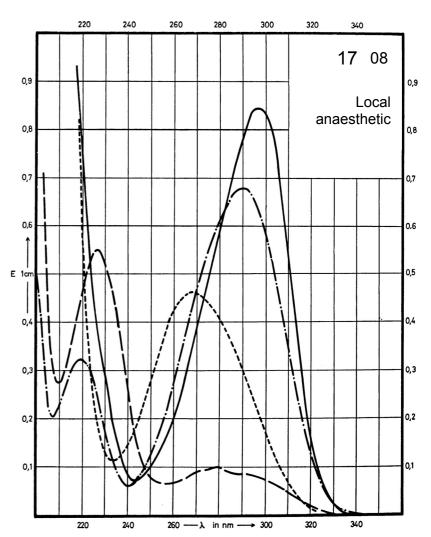

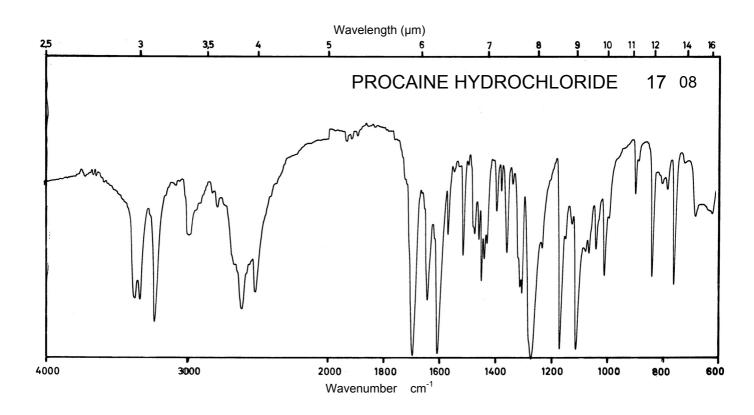





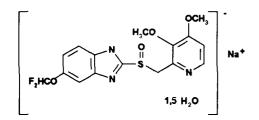


425.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm	291 nm	Decom- position observed	295 nm
E ^{1%} 1cm	391	346		418
ε	16600	14700		17700

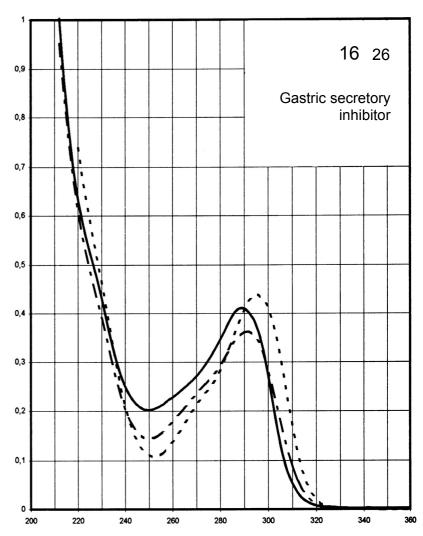

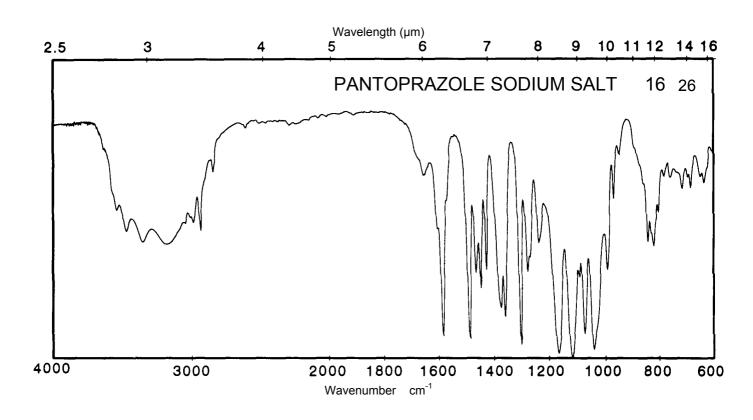





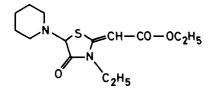


 M_r 298.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 245 nm		281 nm	Decom- position observed
E ^{1%} 1cm	700 296		697	
3	20900 8820		20800	

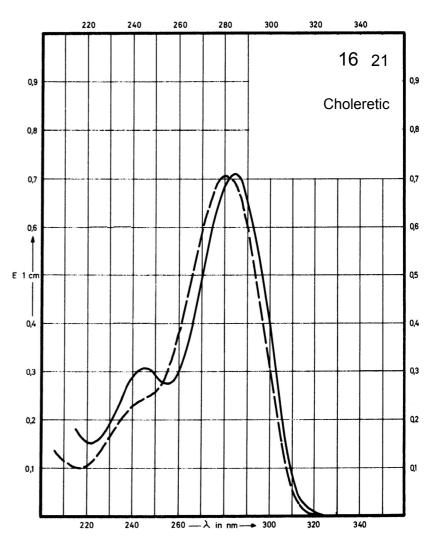

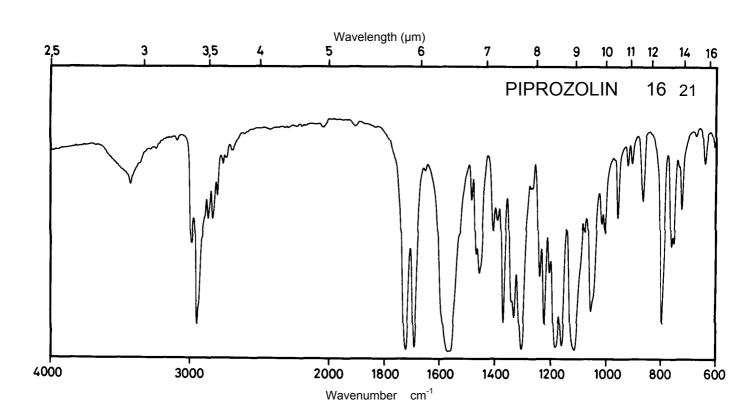



Name PIRENZEPINE DIHYDROCHLORIDE

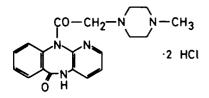


 M_r 424.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm		280 nm	297 nm
E ^{1%} 1cm	186		197	233
ε	7900		8360	9890

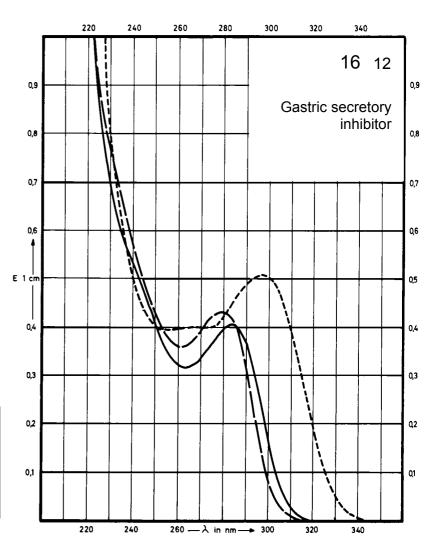

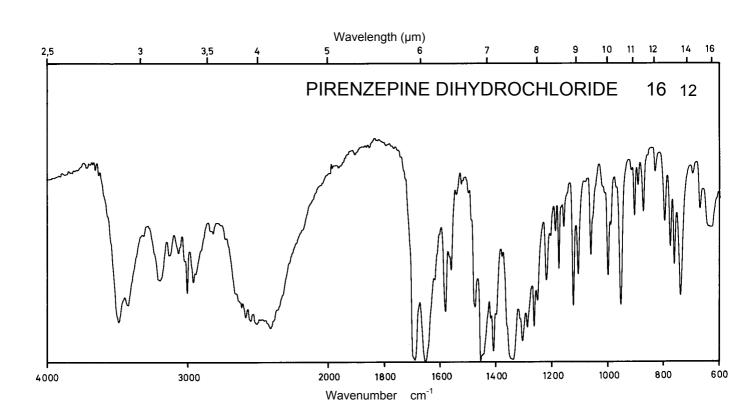


Name PHENOLPHTHALEIN

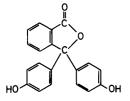


 M_{r}

318.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 228 nm		274 nm 228 nm	294 nm 243 nm
E ^{1%} 1cm	155 938		142 931	175 773
3	4930 29860		4520 29630	5570 24600

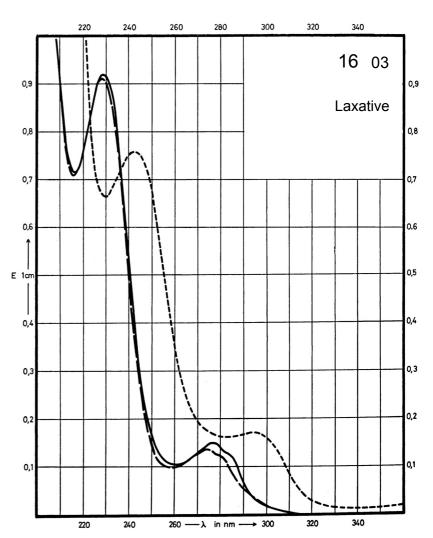

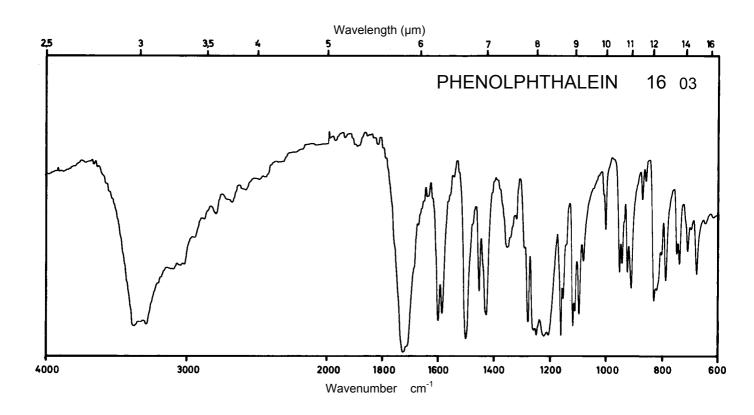







218.3

Concentration


60 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 251 nm			264 nm 258 nm 251 nm
E ^{1%} 1cm	8.6 11.0 9.8			9.2 12.0 11.4
ε	188 240 214			201 262 249

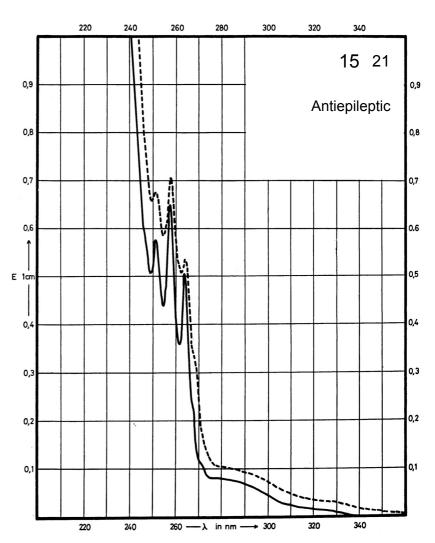

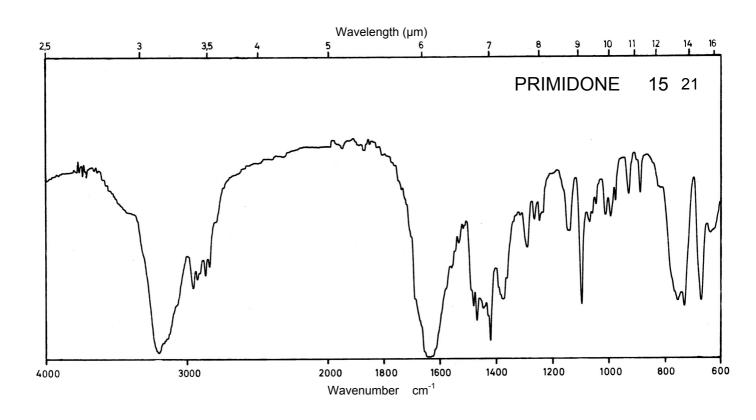





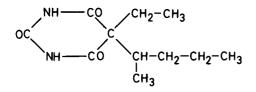


 M_r 226.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				243 nm
E ^{1%} _{1cm}				338
3				7640

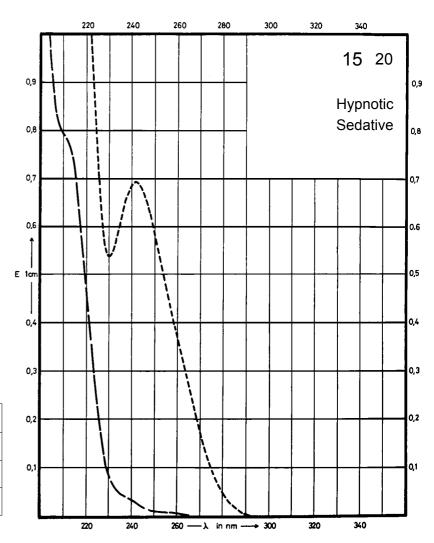

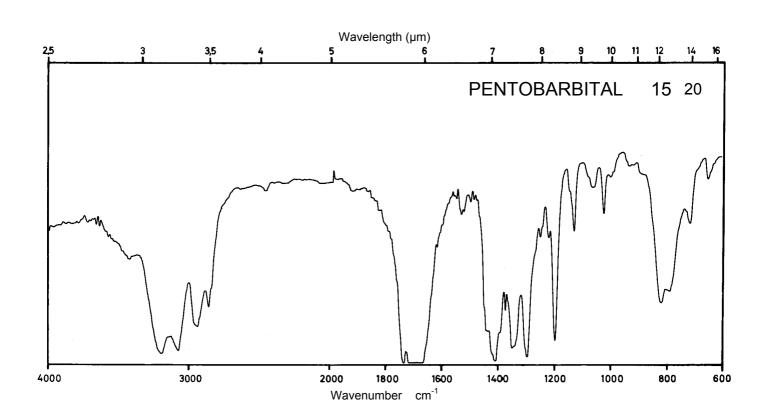



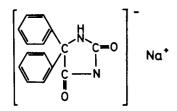



 M_r

Concentration 10 mg / 100 ml


274.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm		259 nm	
E ^{1%} 1cm	15.9		26.5	
ε	435		727	

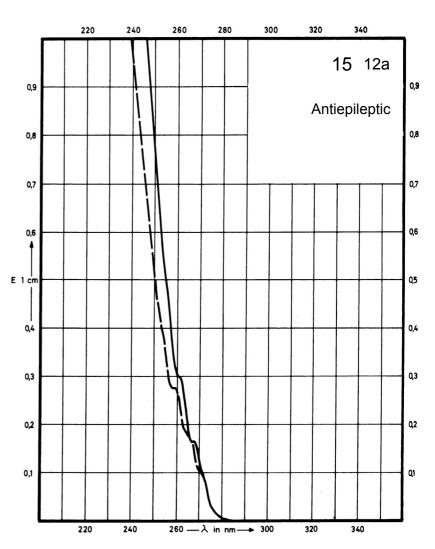

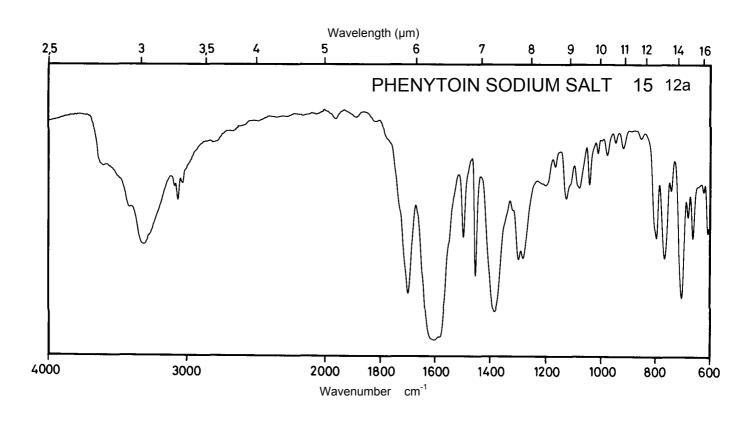

Name PHENYTOIN

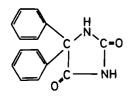


 M_{r}


252.3

Concentration


10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm			
E ^{1%} 1cm	16.0 26.8			
ε	400 680			

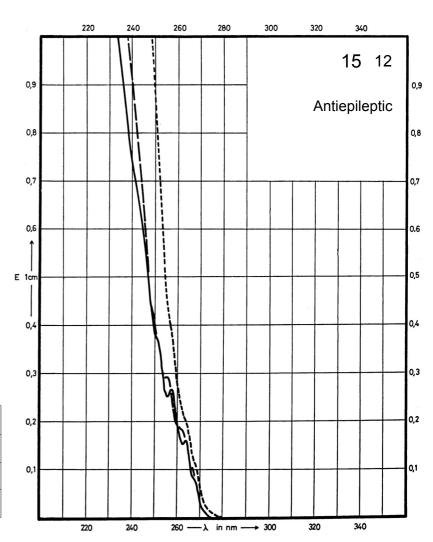

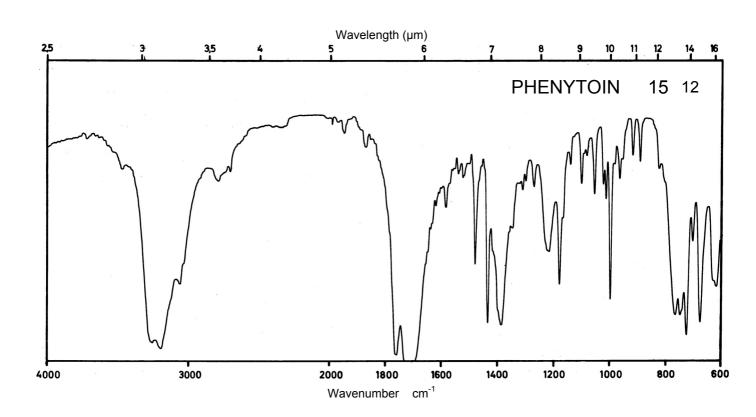





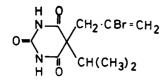


289.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				242 nm
E ^{1%} 1cm				205
ε				5930

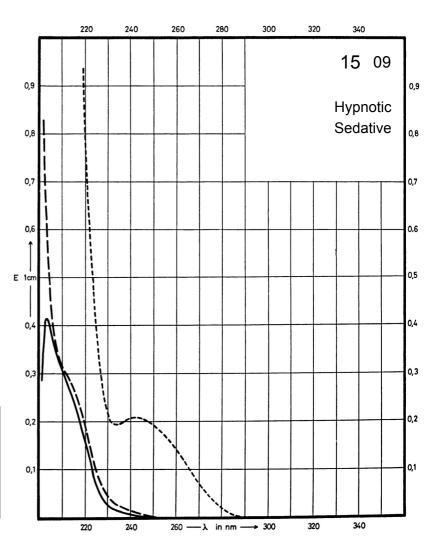

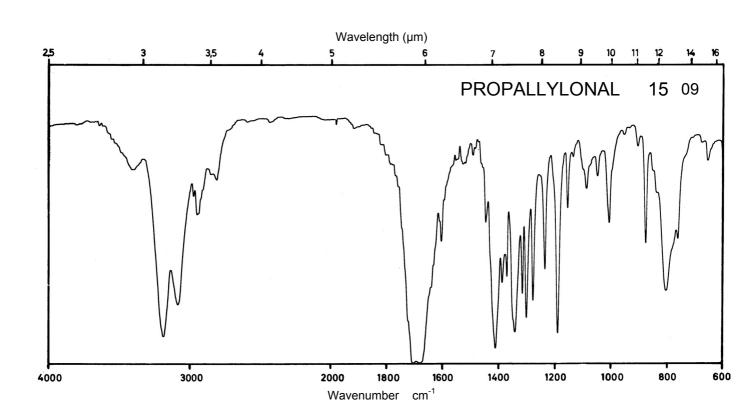

Name PHENOBARBITAL

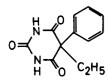



 M_{r}

232.2


Concentration 2 mg / 100 ml 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm			252 nm
E ^{1%} 1cm	26			325
3	600			7550

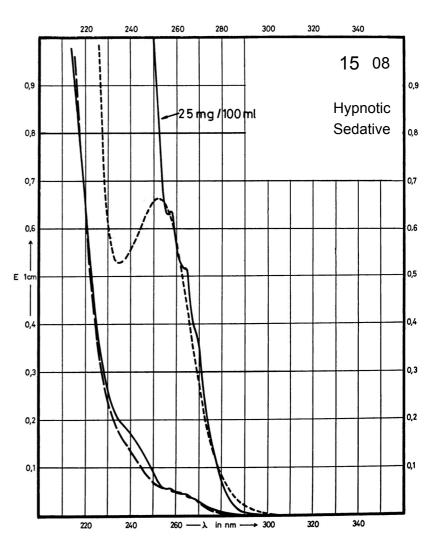

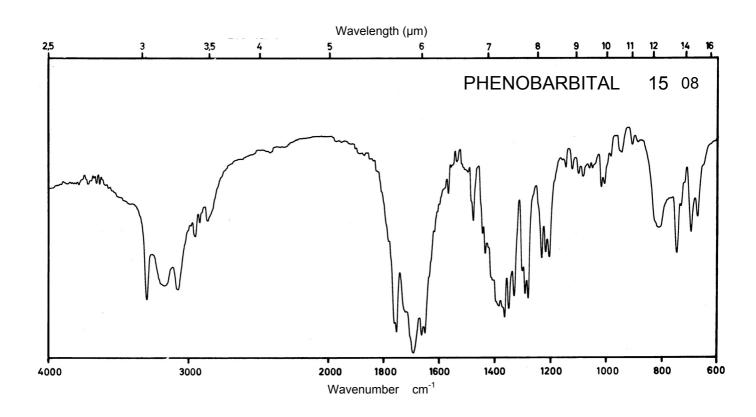





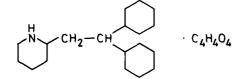


393.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

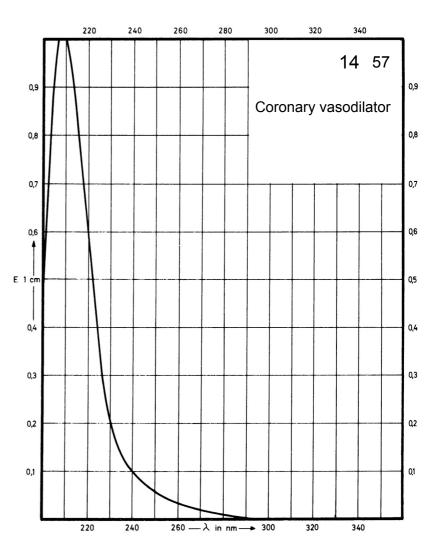

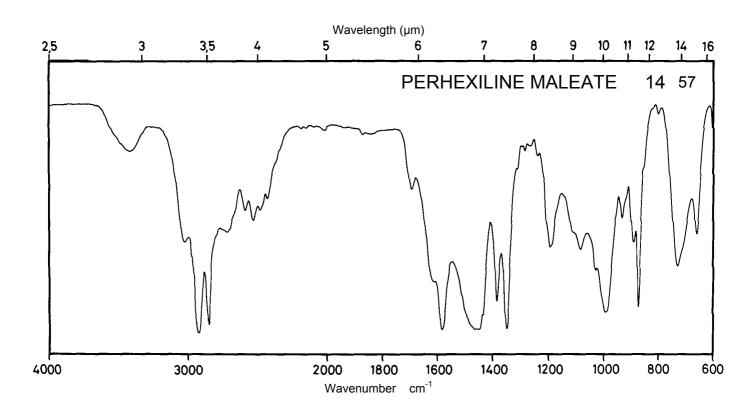





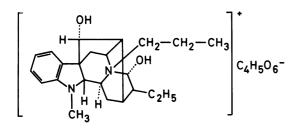


518.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 244 nm		290 nm 245 nm	290 nm 246 nm
E ^{1%} 1cm	53 168		53 160	52 158
ε	2750 8700		2770 8310	2680 8210

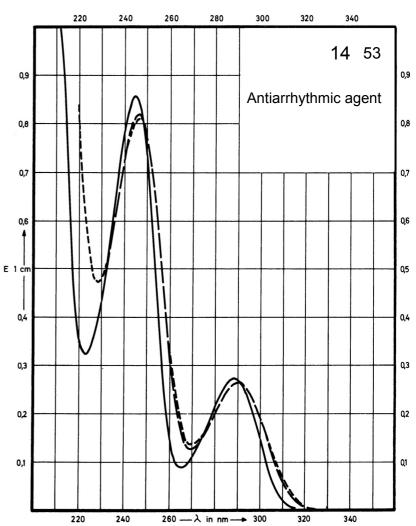

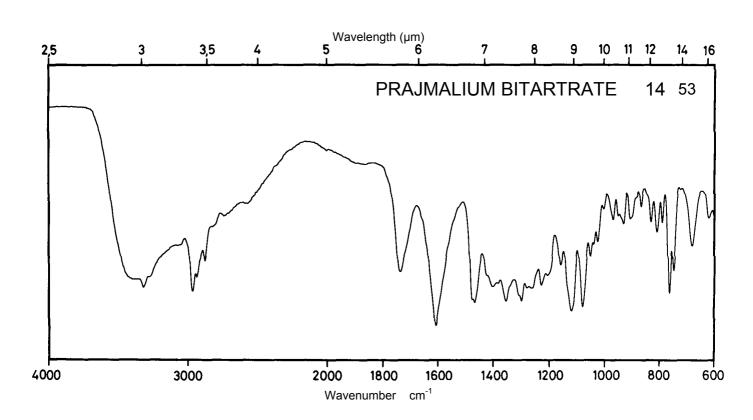



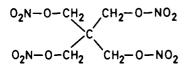



316.1

Concentration


ca. 30 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	256 nm			
E ^{1%} 1cm	ca. 13			
3	ca. 400			

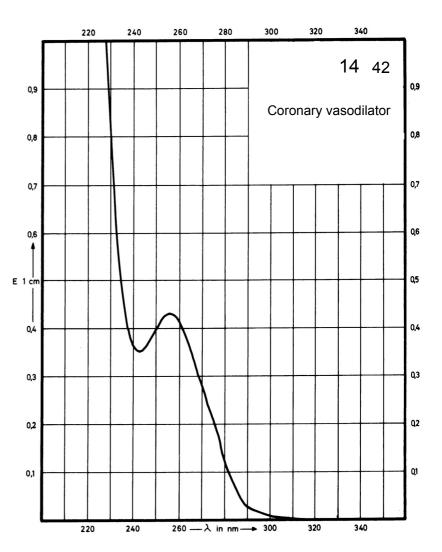

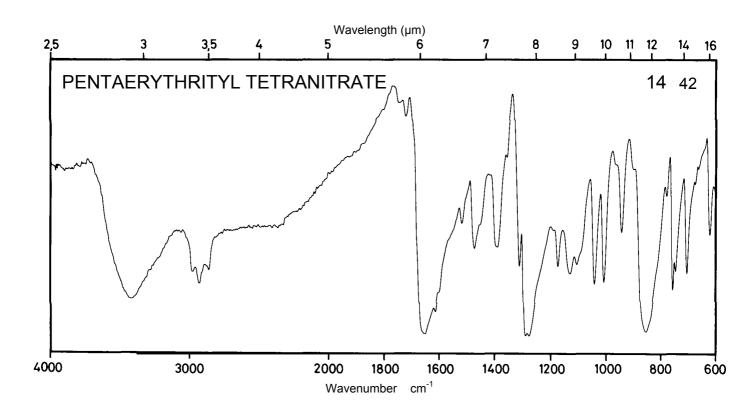





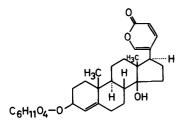

530.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	300 nm	301 nm		Decom- position observed
E ^{1%} 1cm	108	120		
ε	5750	6390		

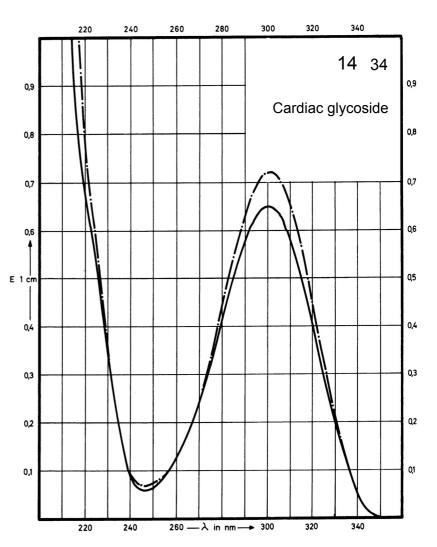

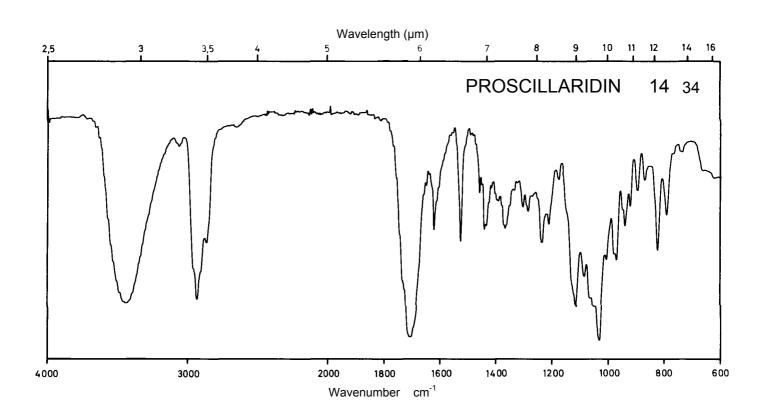



Name PROPAFENONE HYDROCHLORIDE

 M_{r}


377.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	304 nm 248 nm		305 nm 251 nm	306 nm 251 nm
E ^{1%} 1cm	84 217		81 205	81 184
ε	3170 8180		3060 7740	3060 6960

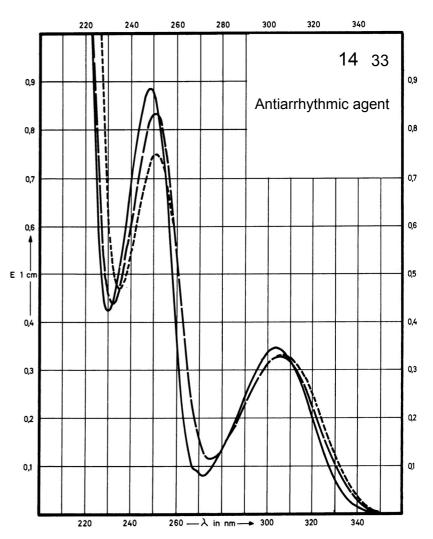

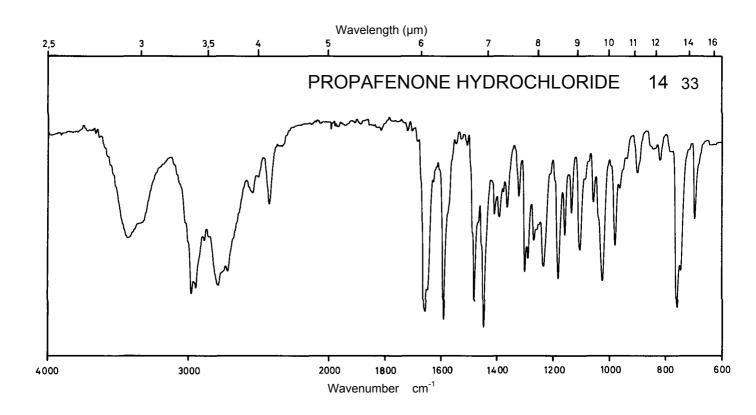





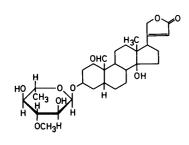

 M_r 548.7

Concentration 2 mg / 100 ml 100 mg / 100 ml

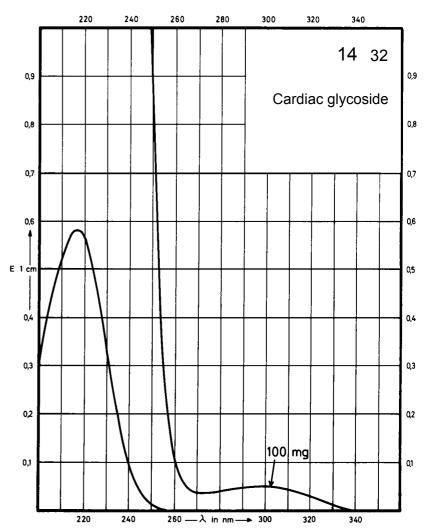

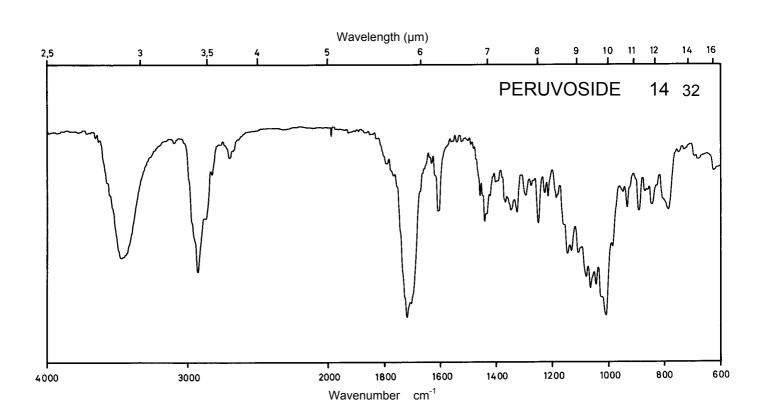


Name PROCAINAMIDE HYDROCHLORIDE


$$\begin{bmatrix} C_{2}H_{5} \\ I \\ I \\ I \\ I \\ C_{2}H_{5} \end{bmatrix}^{+} CI^{-}$$

 M_{r}

271.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm		222 nm	273 nm
E ^{1%} 1cm	680		410	595
ε	18480		11140	16170

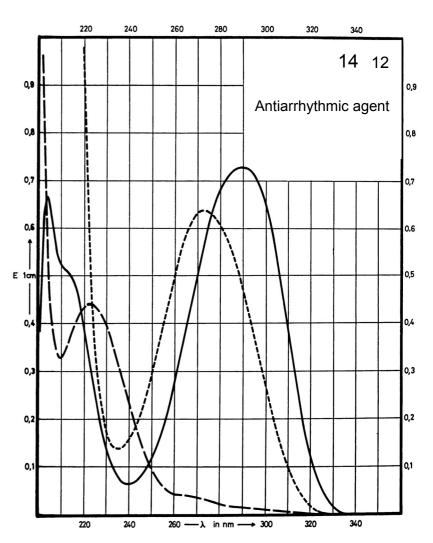


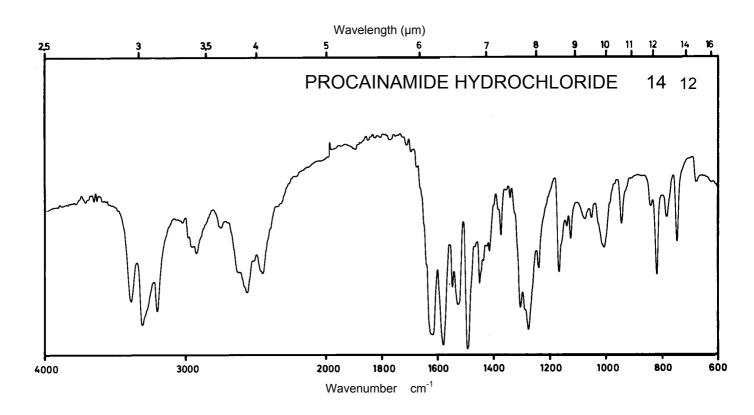





419.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	258 nm	258 nm		
E ^{1%} 1cm	15.0	15.3		
3	630	640		

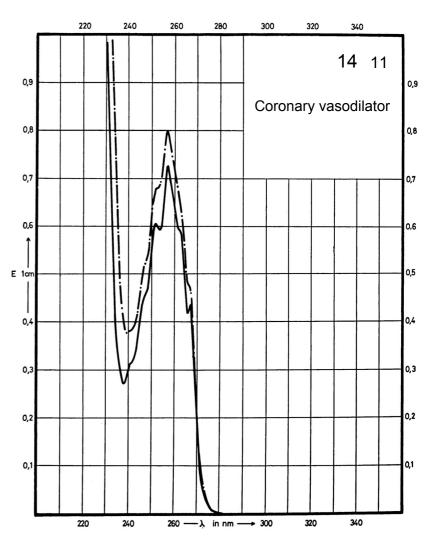

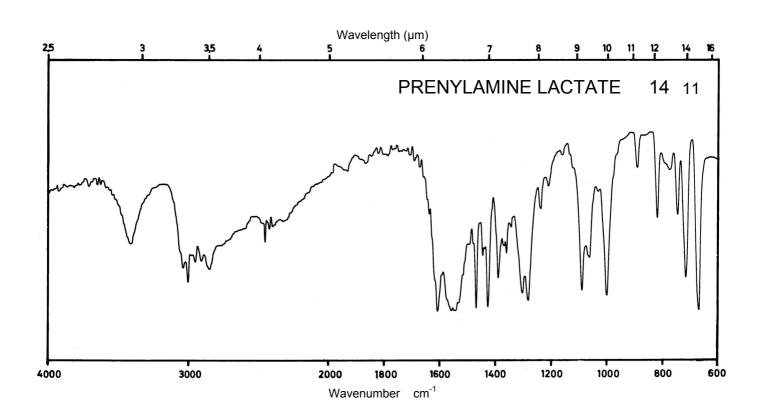





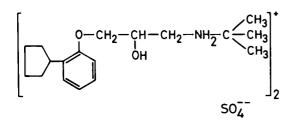

681.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm	269 nm	269 nm	
E ^{1%} 1cm	54	51	51	
3	3680	3470	3470	

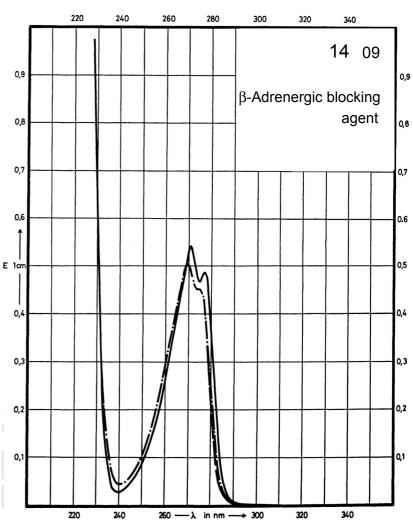

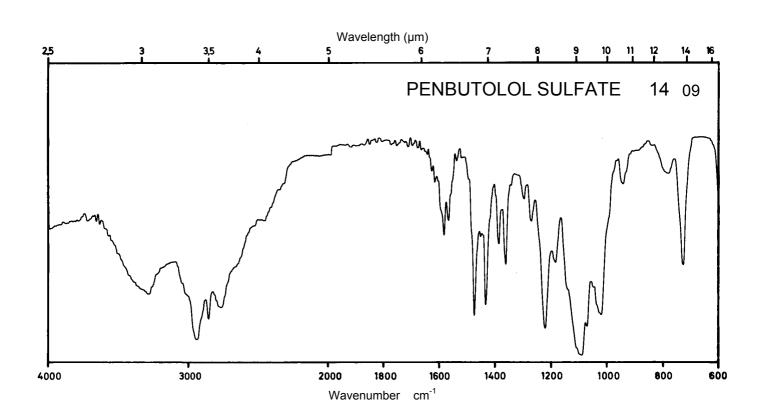


Name POLYTHIAZIDE


$$\begin{array}{c|c} CI & H \\ & CH_2-S-CH_2-CF_3 \\ & CH_3 \end{array}$$

 M_r 439.9

Concentration 1 mg / 100 ml 6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	314 nm 270 nm 226 nm		314 nm 271 nm 227 nm	Decom- position observed
E ^{1%} 1cm	68 497 891		71 484 900	
ε	3000 21860 39200		3120 21290 39560	

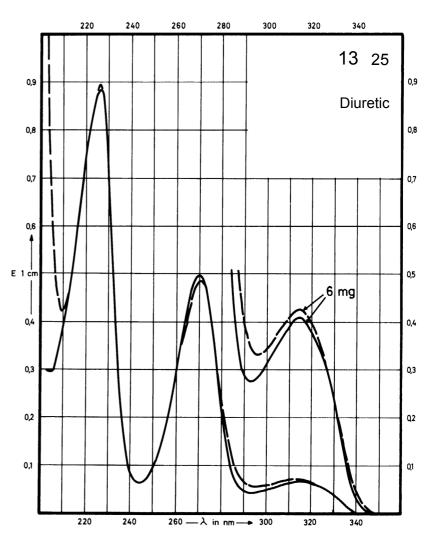

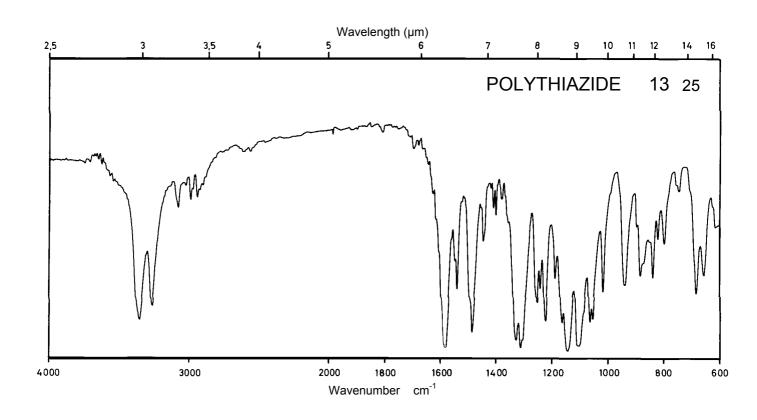

Name PIRETANIDE

H₂NO₂S COOH


 M_{r}


362.4

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	348 nm 275 nm	345 nm 275 nm	355 nm 273 nm 266 nm	330 nm 267 nm
E ^{1%} 1cm	90 313	78 273	31 131 132	76 285
ε	3260 11343	2830 9890	1120 4750 4780	2750 10330

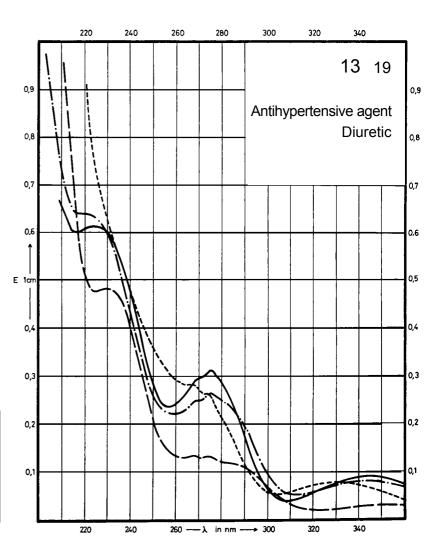

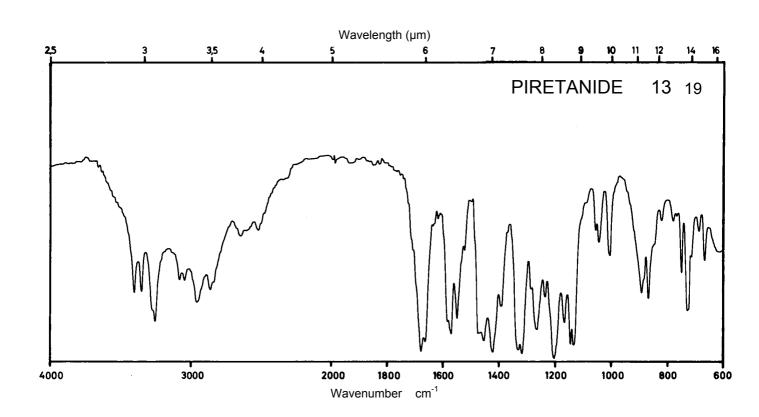





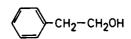

122.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	252 nm 258 nm 268 nm			
E ^{1%} 1cm	13.3 16.0 10.7			
ε	163 196 131			

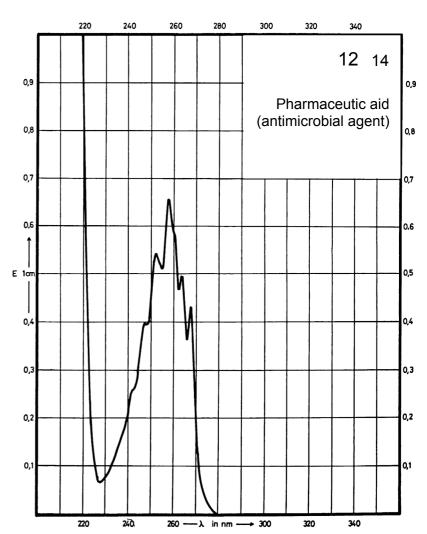

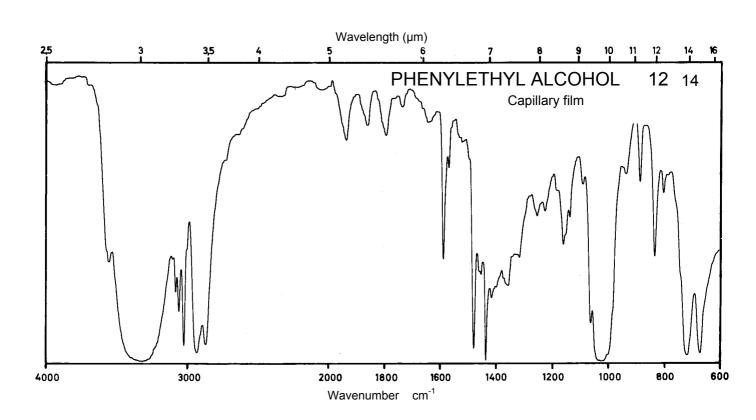


Name PHENOL


OH

 M_{r}

94.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		270 nm	287 nm
E ^{1%} 1cm	201		172	287
ε	1890		1620	2700

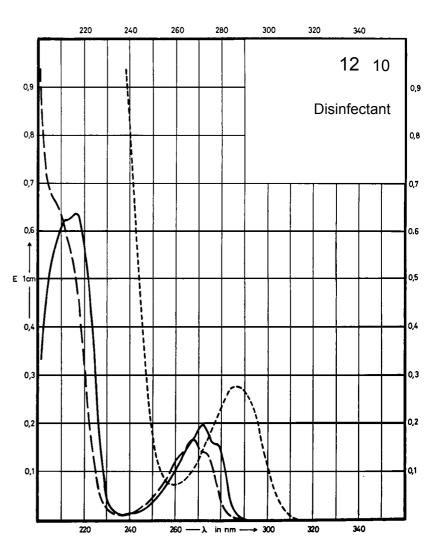

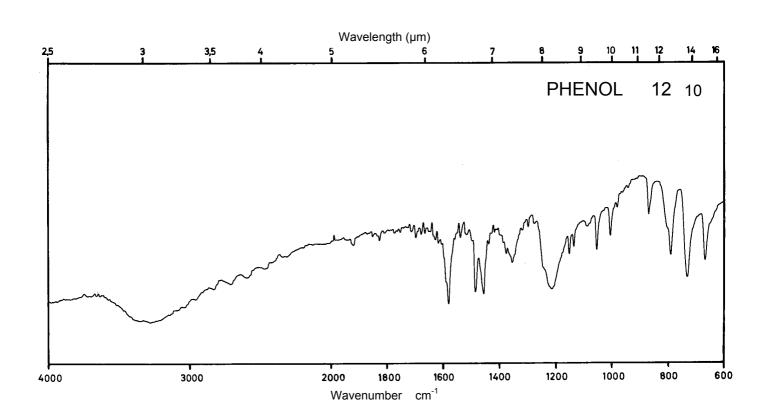

Name PROPYLPARABEN SODIUM SALT


 M_{r}


202.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	256 nm		255 nm	296 nm
E ^{1%} _{1cm}	827		789	1208
ε	16550		15800	24200

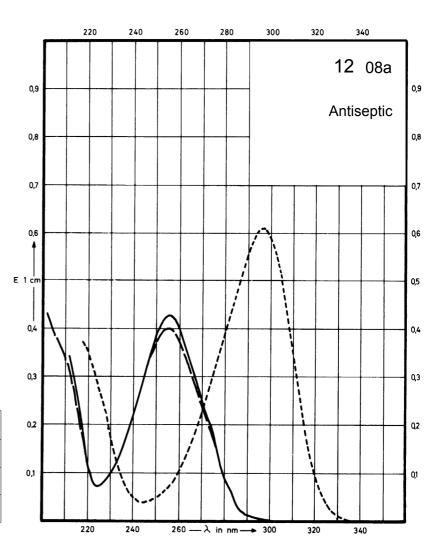

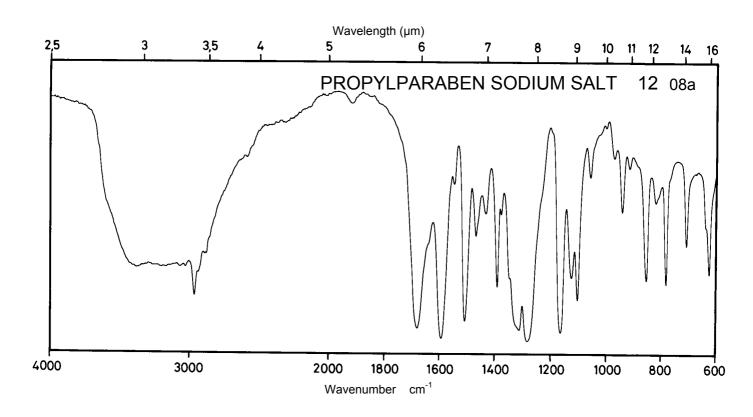



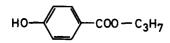



180.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	256 nm		255 nm	296 nm
E ^{1%} 1cm	928		877	1324
ε	16720		15800	23860

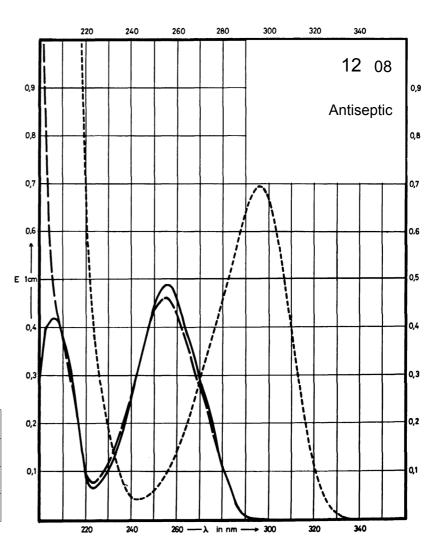

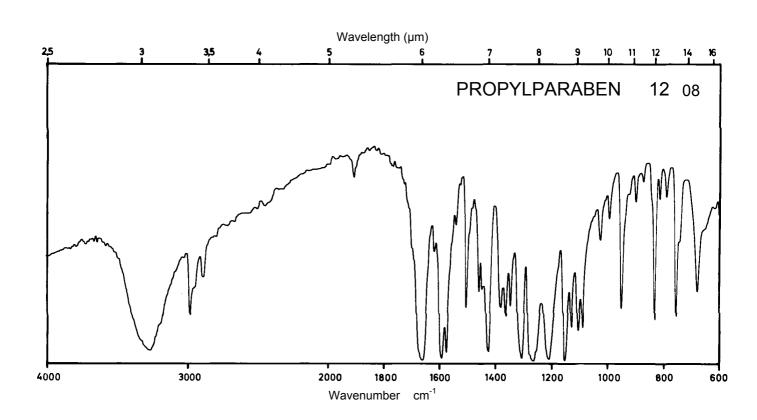



Name PIROCTONE OLAMINE

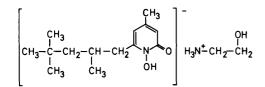


 M_r 298.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	304 nm 229 nm	308 nm	294 nm 224 nm	312 nm
E ^{1%} 1cm	234 240	214	236 143	225
ε	6980 7170	6390	7040 4270	6720

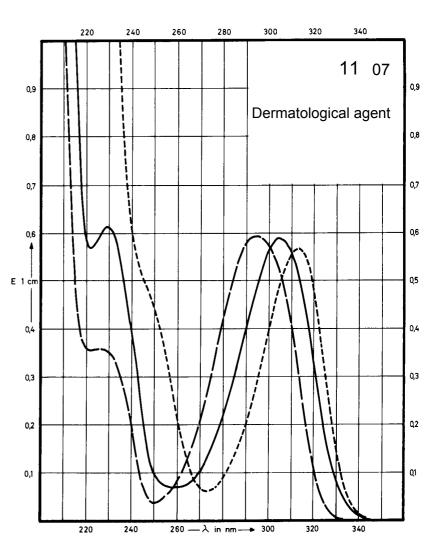

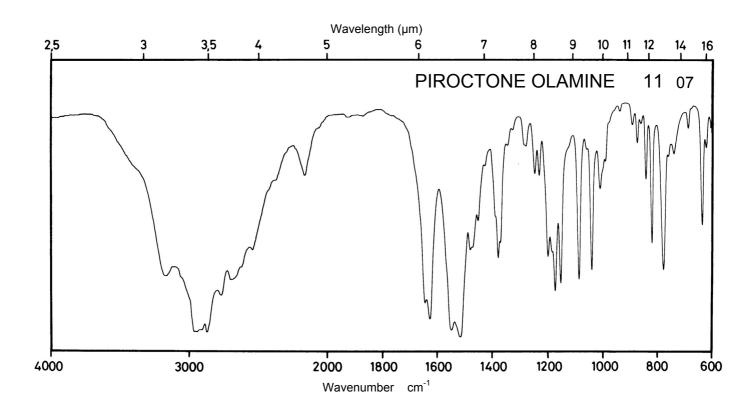





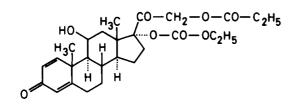


 M_r 488.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	242 nm			
E ^{1%} 1cm	317			
3	15500			

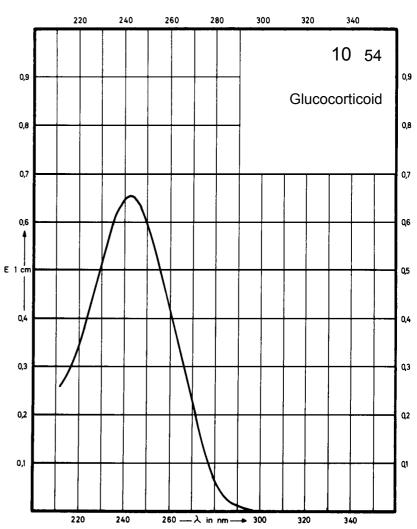

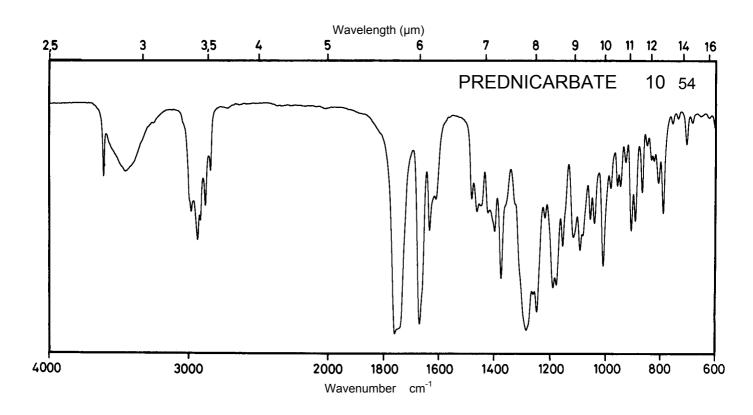


Name PARAMETHASONE 21-ACETATE

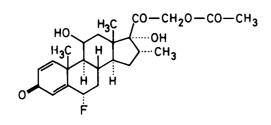


 M_r


434.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm	245 nm		
E ^{1%} 1cm	368	380		
3	16000	16530		

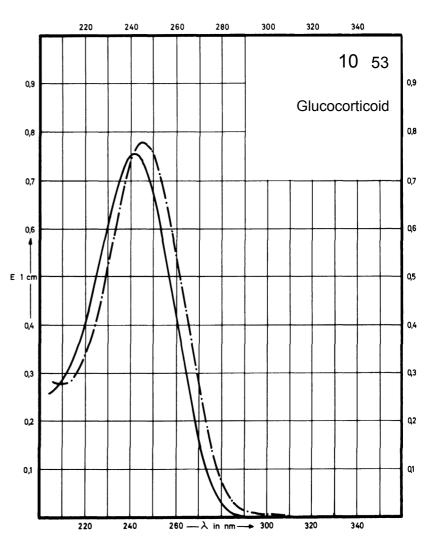

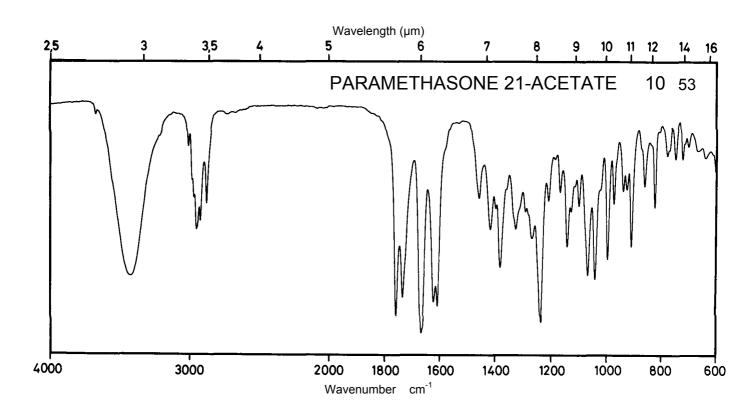



Name PREDNYLIDENE

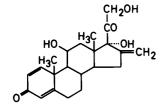


 M_r 372.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	244 nm	248 nm		
E ^{1%} 1cm	412	418		
3	15350	15570		

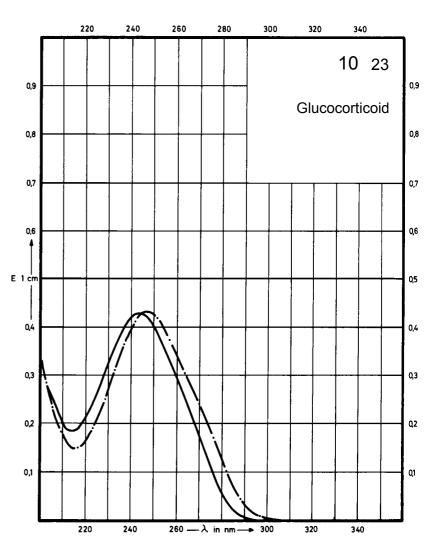

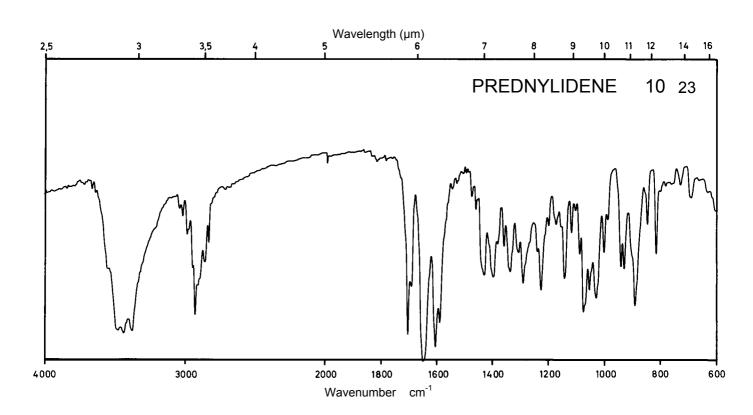



Name PREDNISOLONE 21-ACETATE

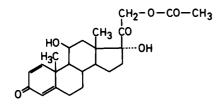


 M_r 402.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm	248 nm		
E ^{1%} 1cm	370	370		
3	14890	14890		

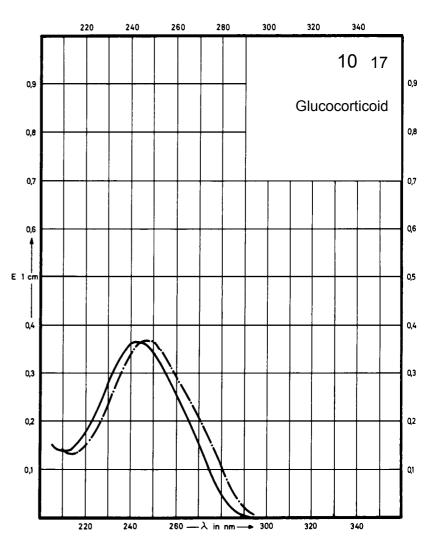

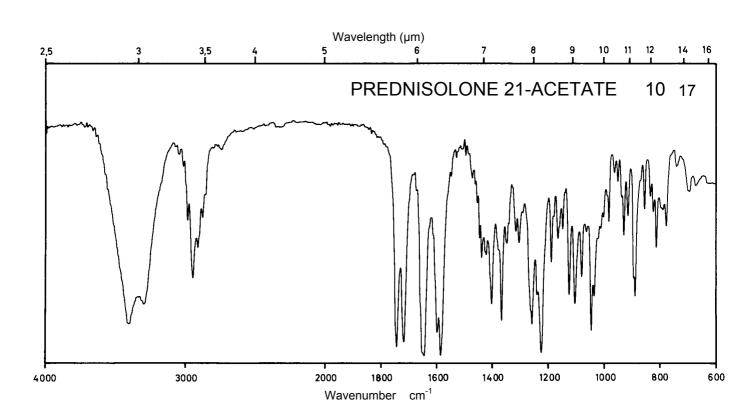





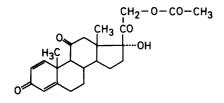


 M_{r} 400.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm	244 nm		
E ^{1%} 1cm	380	390		
3	15220	15620		

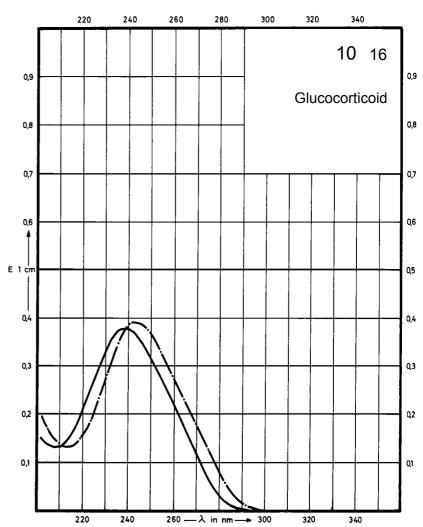

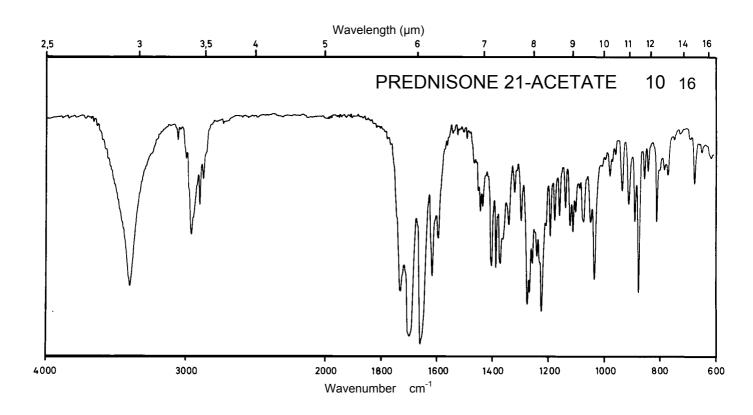




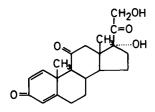


 M_{r}


358.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm			
E ^{1%} 1cm	421			
3	15090			

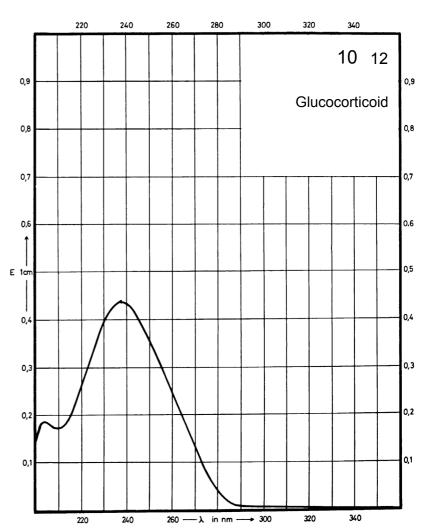

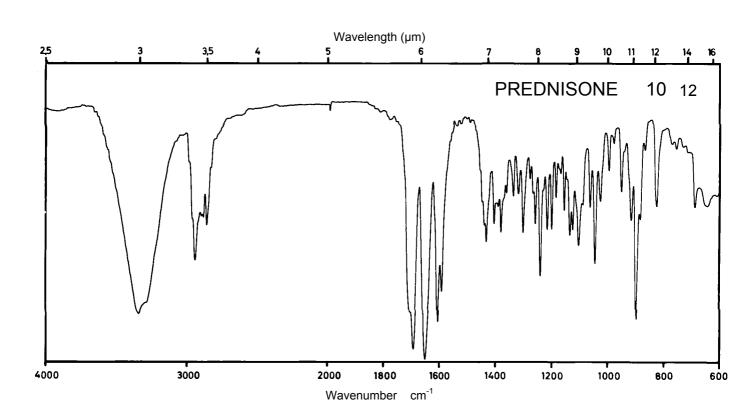





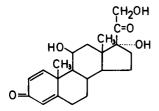


 M_r 360.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	242 nm			
E ^{1%} 1cm	416			
ε	15000			

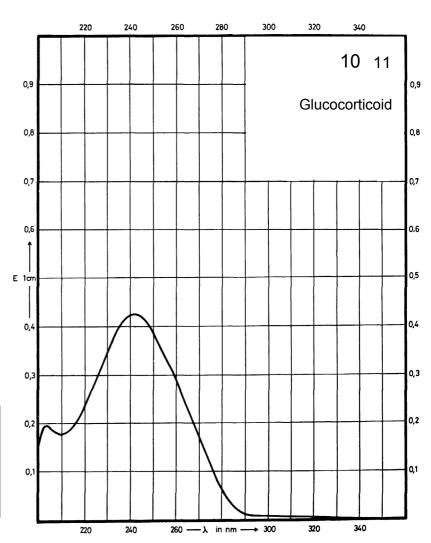

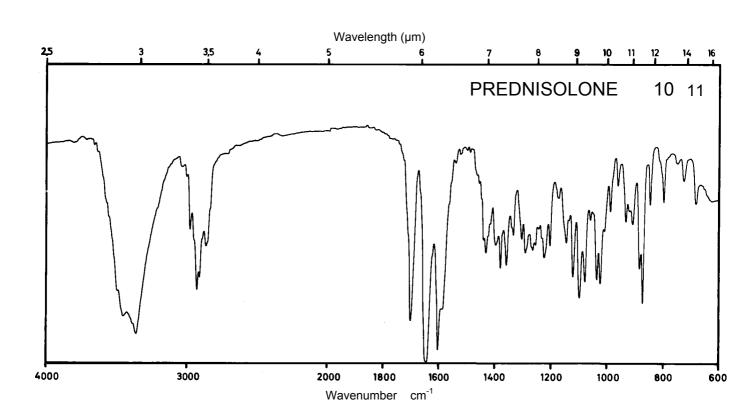



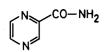



 M_r

123.1


Concentration 0.75 mg / 100 ml 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	320 nm	310 nm	310 nm	310 nm
	268 nm	268 nm	268 nm	268 nm
E ^{1%}	43	50	50	50
1cm	640	658	658	658
3	530	620	620	620
	7900	8100	8100	8100

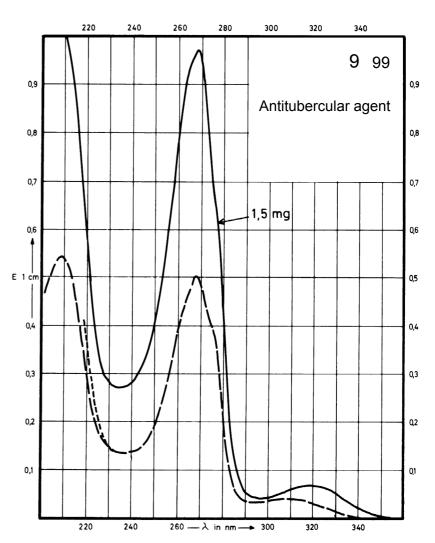

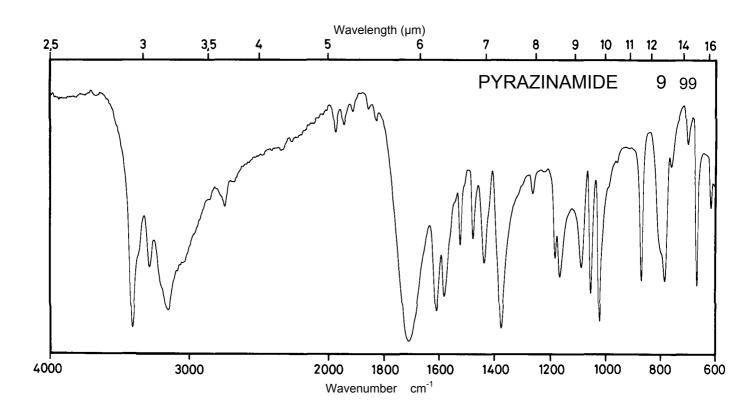



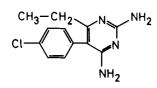



 M_r

Concentration 1 mg / 100 ml


248.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	286 nm	285 nm	274 nm	287 nm
E ^{1%} 1cm	365	340	320	381
3	9080	8450	7960	9480

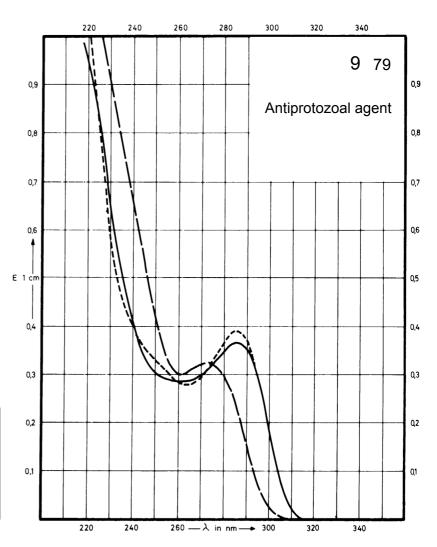

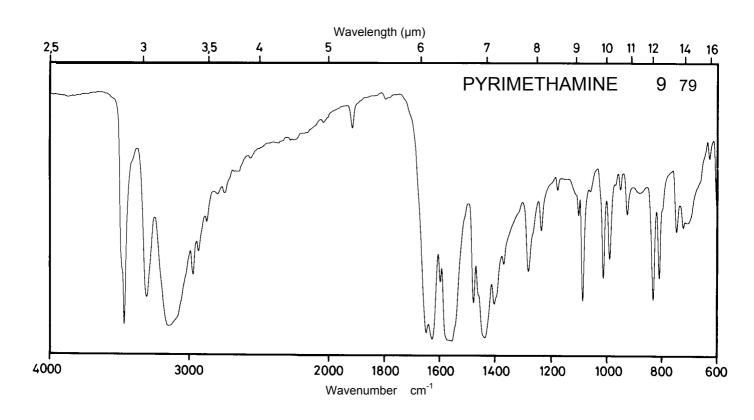



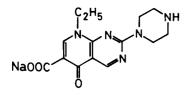



325.3

Concentration


0.3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	325 nm 276 nm		321 nm 275 nm	332 nm 272 nm
E ^{1%} 1cm	300 1795		340 2150	400 1720
ε	9800 58400		11000 69800	13000 56000

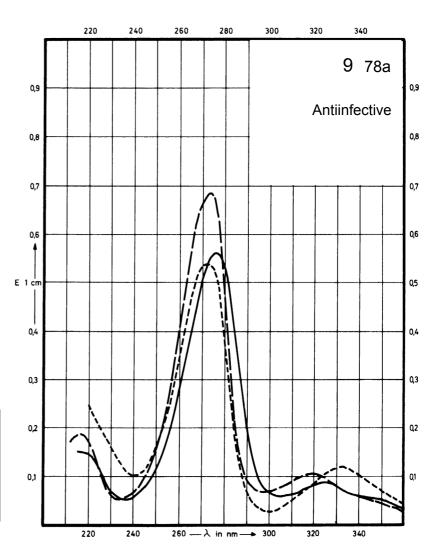

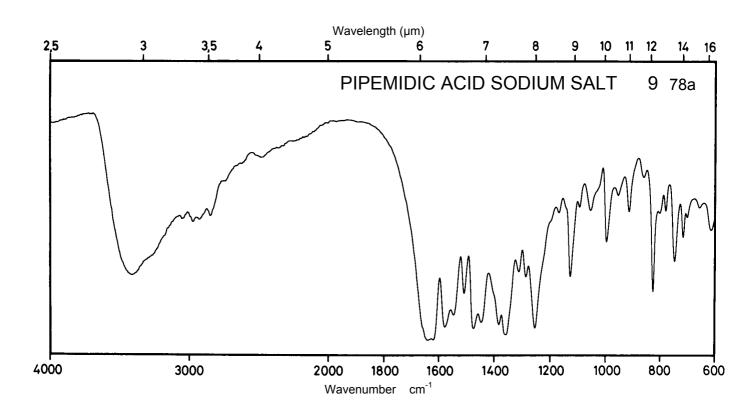



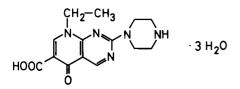



357.4

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	325 nm	330 nm	321 nm	332 nm
	276 nm	264 nm	275 nm	272 nm
E ^{1%}	280	392	306	368
1cm	1640	1525	1930	1550
3	10000	14010	10940	13150
	58600	54500	69000	55400

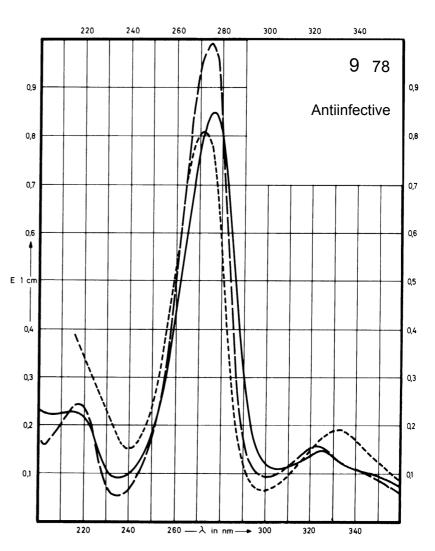

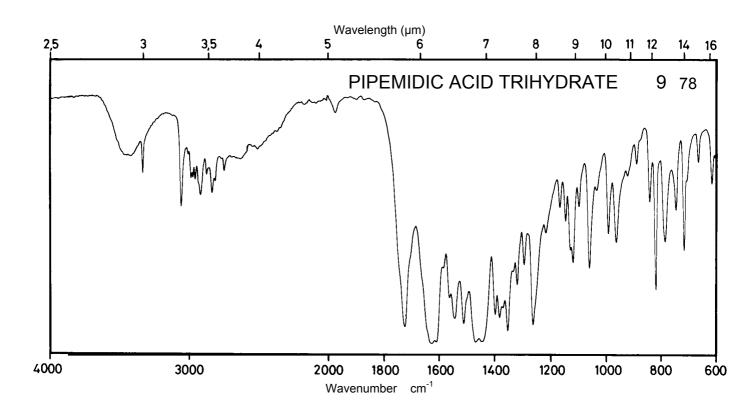





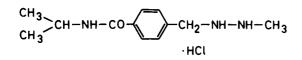


 M_r 257.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	231 nm	232 nm	233 nm	235 nm
E ^{1%} 1cm	Decom- position observed	488	488	500
ε		12580	12580	12890

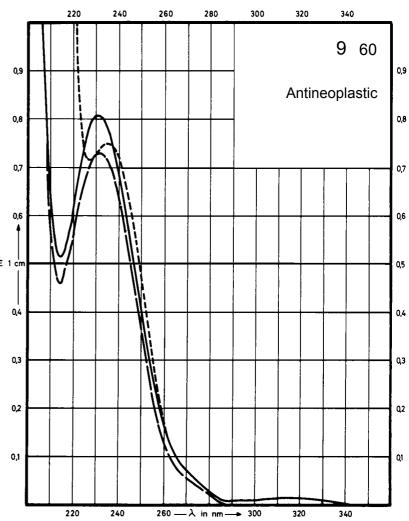

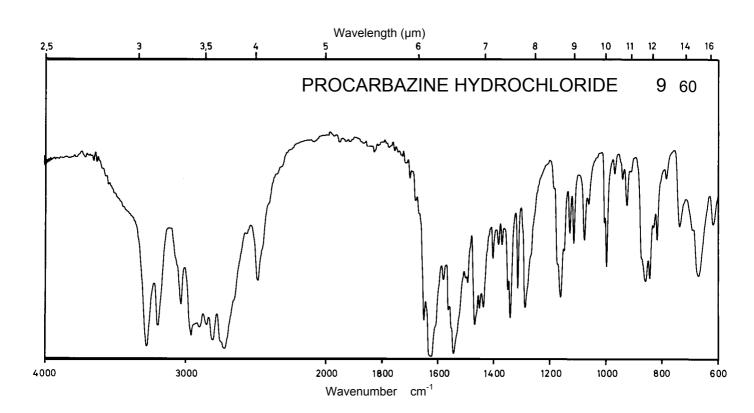



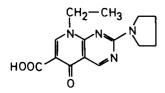



 M_r 288.3

Concentration 0.5 mg / 100 ml 2 mg / 100 ml


Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of absorption			350 nm* 324 nm*	332 nm 274 nm
'			282 nm*	
E ^{1%} 1cm			183* 250*	363
1cm			1720*	1680
			5270*	10480
3			7190* 49600*	48430

^{* 1} M HCI + methanol (1+9)

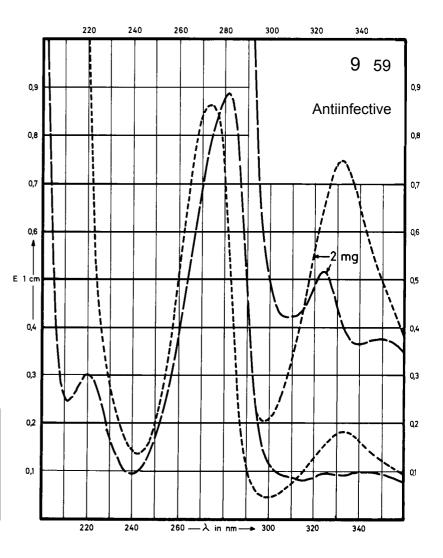

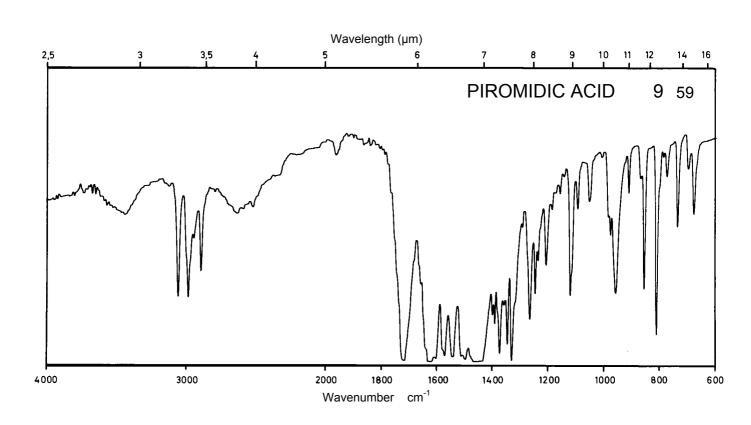





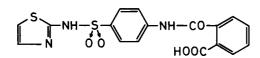



403.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 259 nm		284 nm 250 nm	260 nm
E ^{1%} 1cm	651 503		574 480	585
3	26260 20290		23160 19360	23600

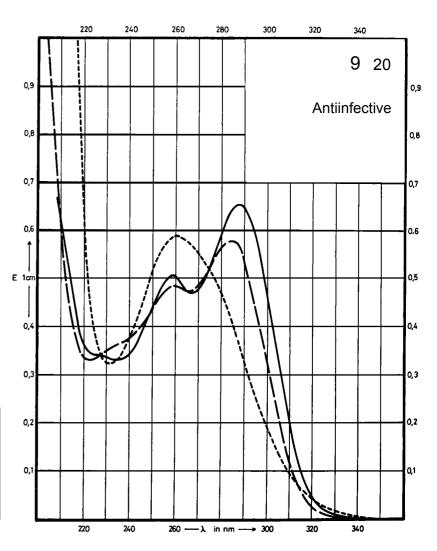


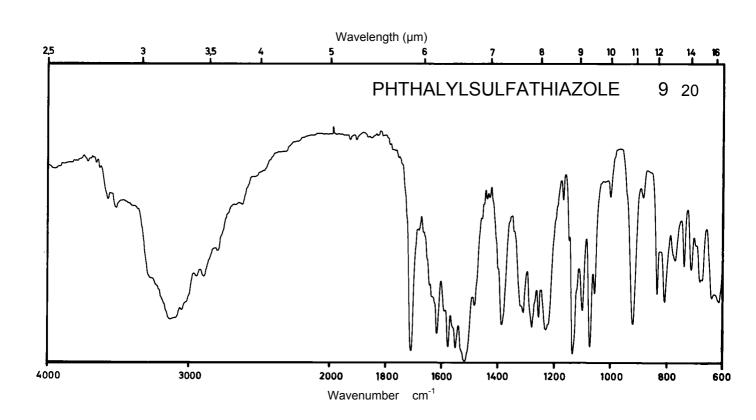

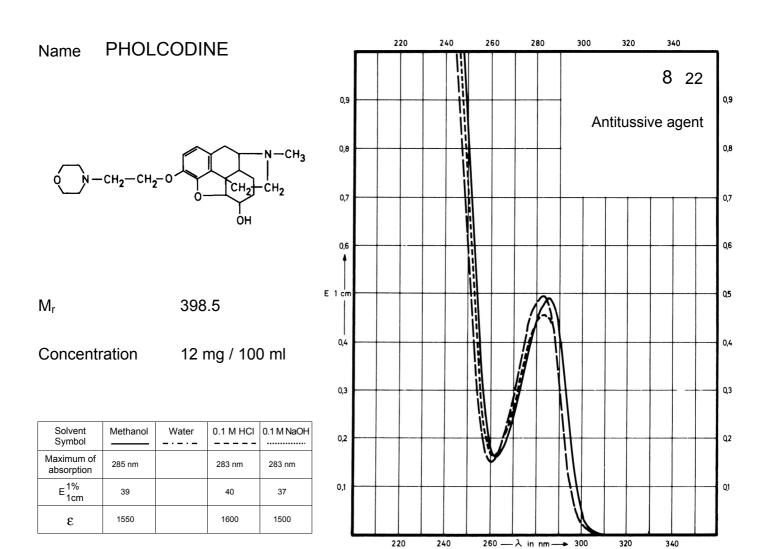

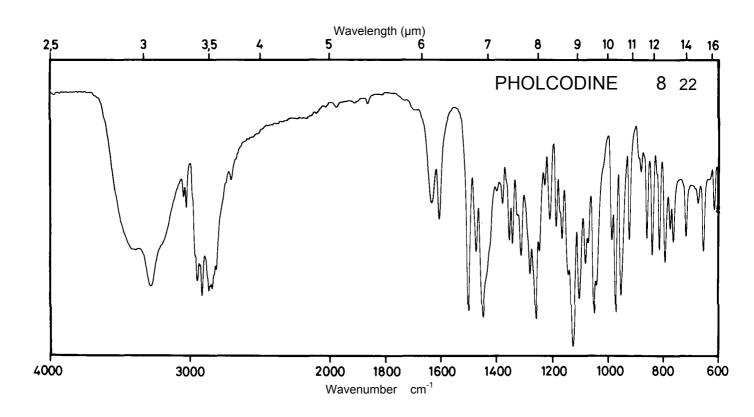





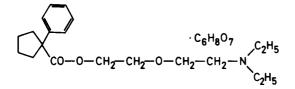

 M_r 525.6

Concentration 10 mg / 100 ml 100 mg / 100 ml

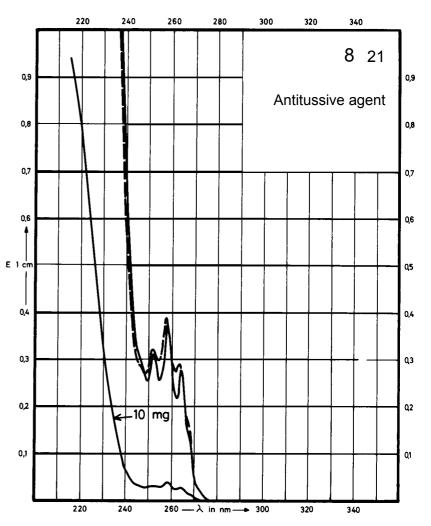

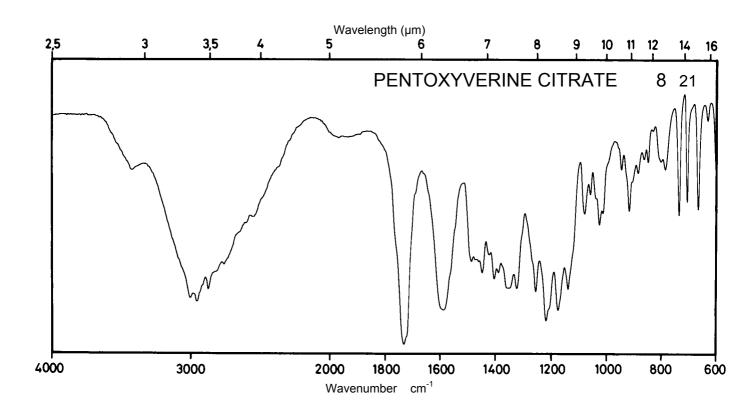



Name PIPAZETATE HYDROCHLORIDE

 M_r 436.0

Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol 274 nm 254 nm 226 nm Maximum of 271 nm 251 nm 271 nm absorption 114 170 E^{1%} 1cm 123 188 122 ε 7420 5370 8190 5330

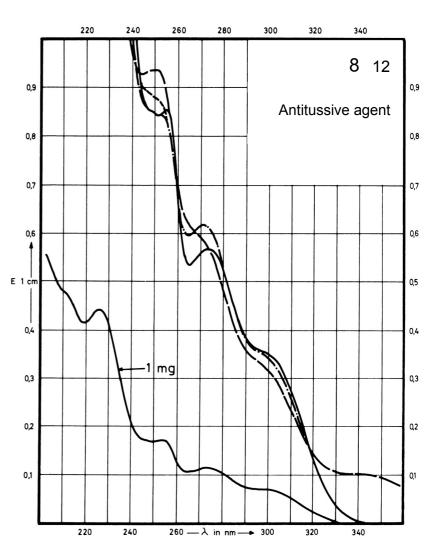

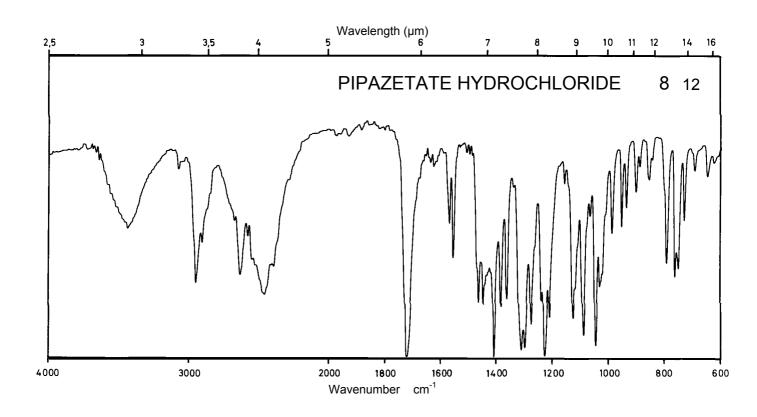





 M_r 1076.6


Concentration 1.2 mg / 100 ml 12 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	318 nm		320 nm	
E ^{1%} 1cm	57		62	
ε	6140		6700	

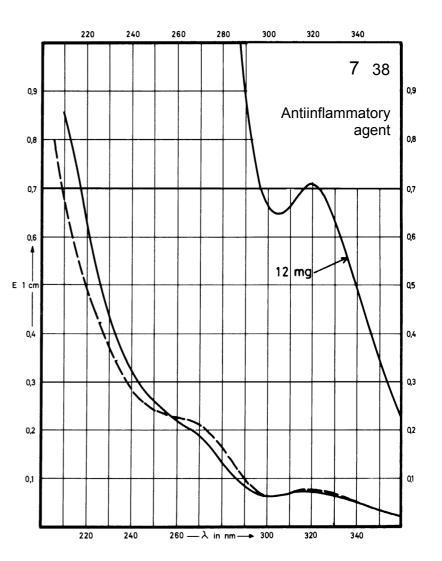

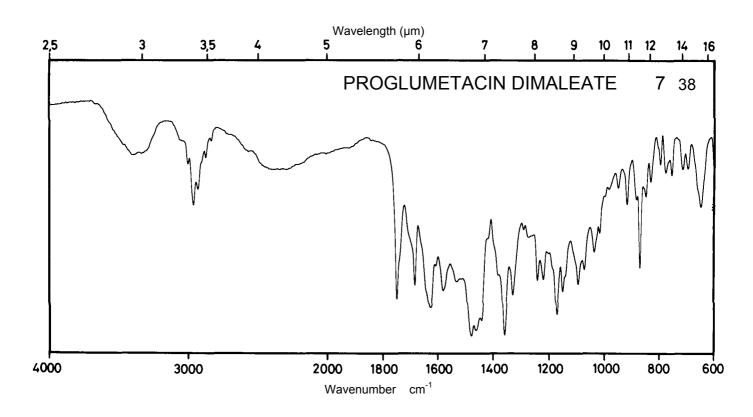





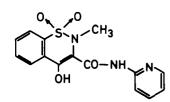

331.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	325 nm 240 nm		334 nm 241 nm	353 nm 286 nm 251 nm
E ^{1%} 1cm	556 368		813 361	475 300 395
ε	18430 12200		26940 11960	15740 9940 13090

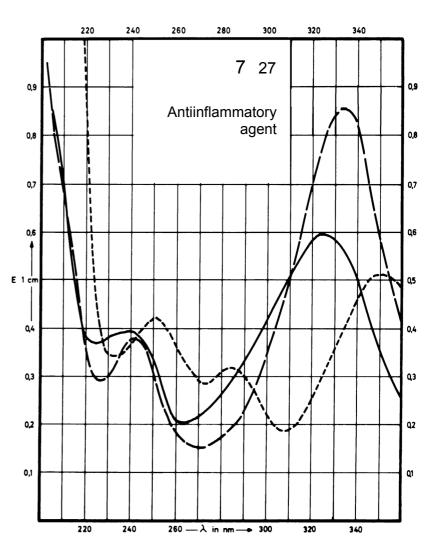


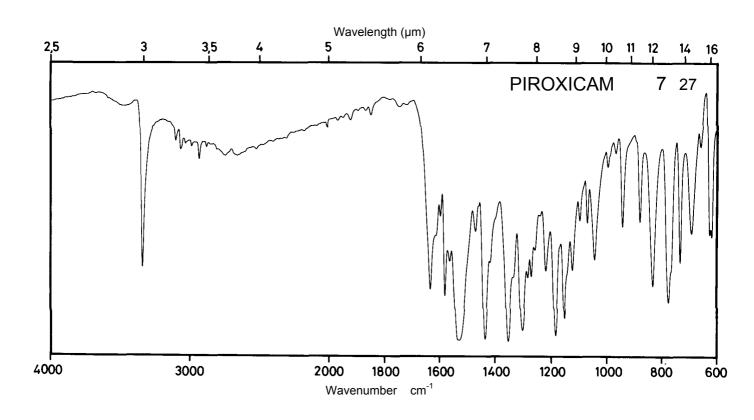


Name PROQUAZONE

 M_r 278.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	356 nm	353 nm	393 nm	353 nm
	276 nm	285 nm	319 nm	285 nm
	235 nm	234 nm	232 nm	234 nm
E ^{1%} 1cm	207 366 1350	218 382 1260	240 532 1155	218 382 1275
3	5750	6080	6680	6080
	10180	10650	14800	10650
	37580	35080	32150	35480

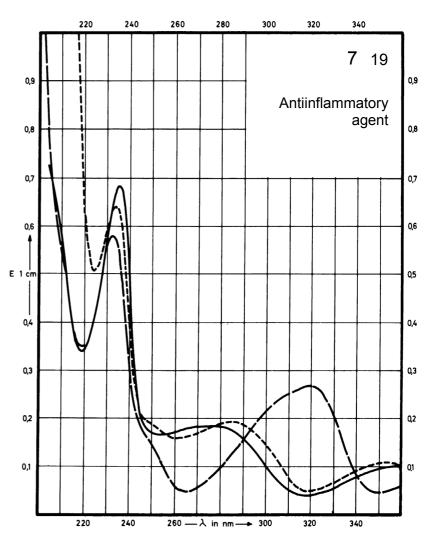

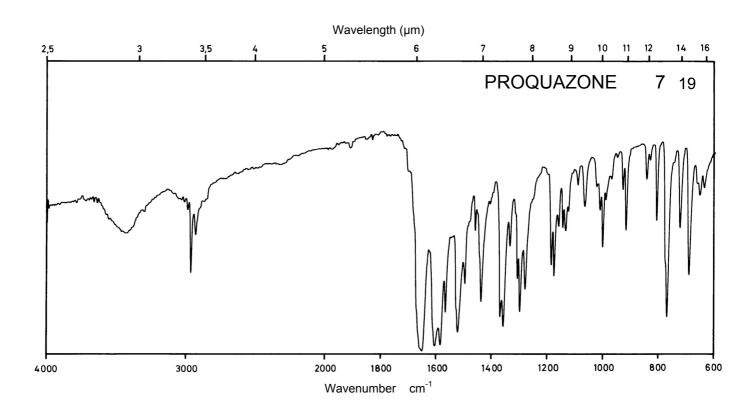

Name PHENYLBUTAZONE


 M_{r}


308.4

Concentration


0.9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm		235 nm	263 nm
E ^{1%} 1cm	482		440	669
ε	14860		13570	20630

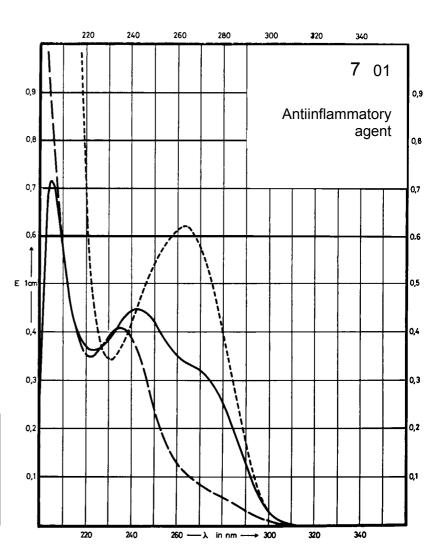

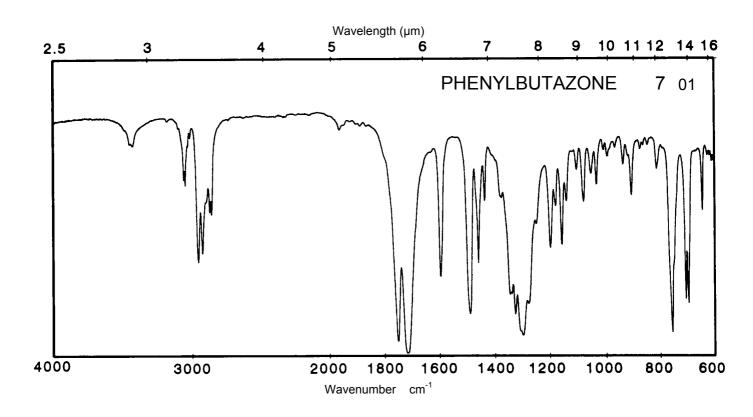




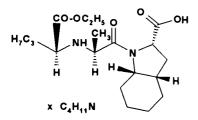


441.6


Concentration 0.5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

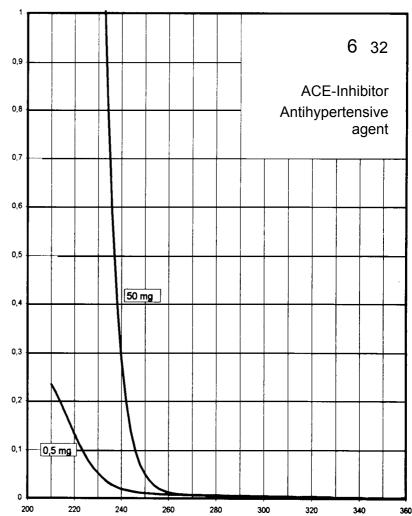

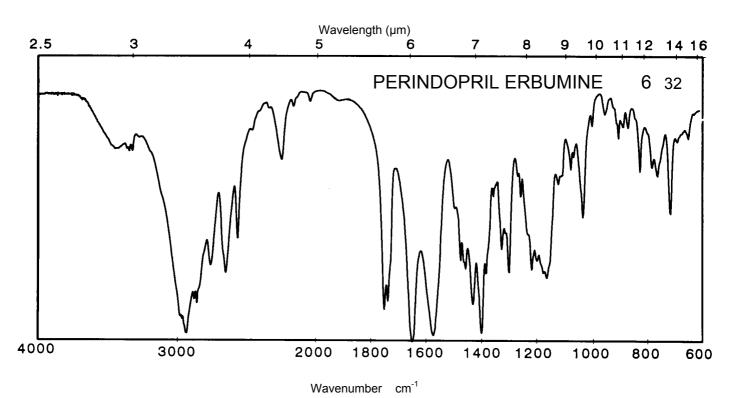


Name PHENOXYBENZAMINE HYDROCHLORIDE

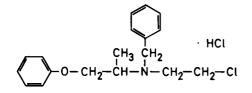


 M_r


340.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 270 nm 260 nm		275 nm 269 nm 263 nm	
E ^{1%} 1cm	34 45 37		31 43 36	
ε	1150 1550 1250		1060 1480 1230	

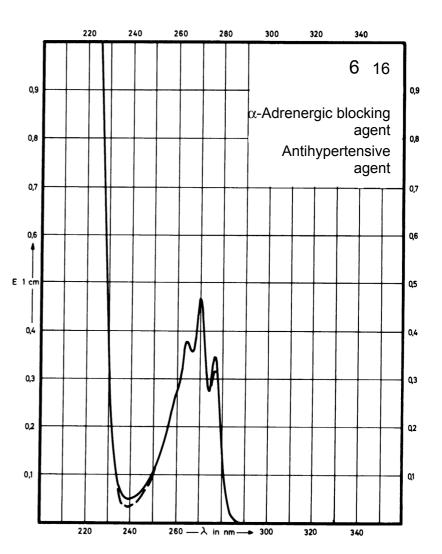

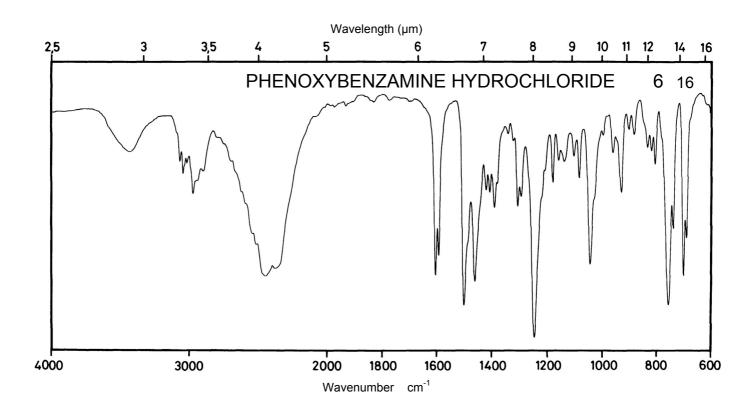

Name PRAZOSIN HYDROCHLORIDE

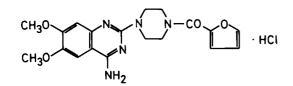



 M_{r}

419.9


Concentration 0.5 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	342 nm 331 nm 250 nm		331 nm 247 nm	345 nm 252 nm
E ^{1%} 1cm	222 233 1430		257 1343	137 1500
ε	9300 9760 60100		10800 56400	5760 62900

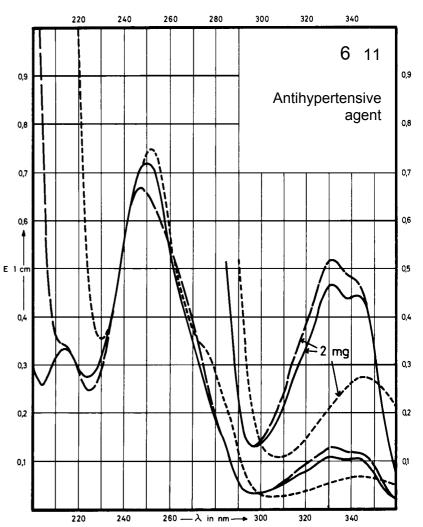

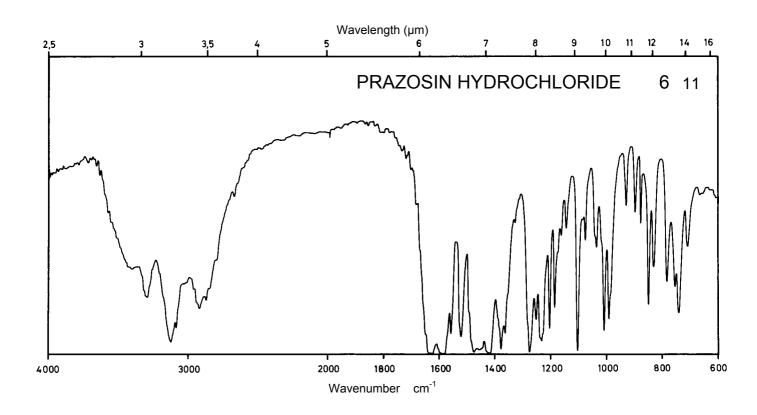

Name PHENTOLAMINE HYDROCHLORIDE

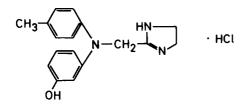


 M_{r}


317.8


Concentration


2.5 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm		278 nm	291 nm
E ^{1%} 1cm	206		226	241
ε	6540		7180	7670

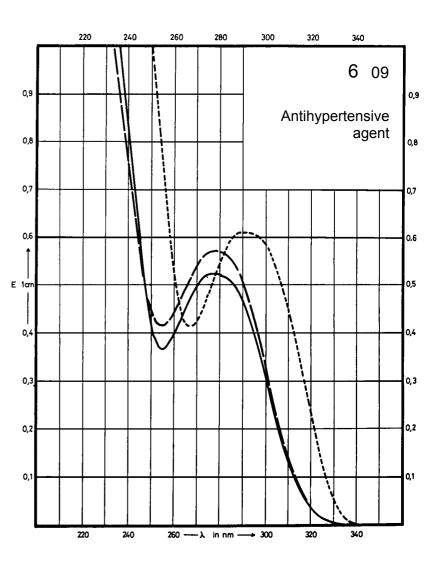

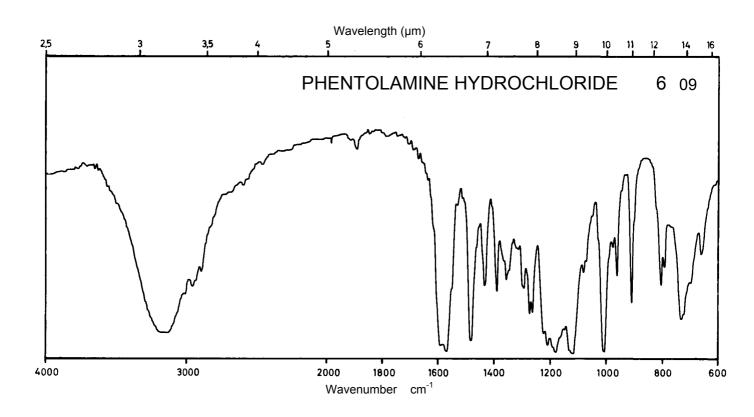

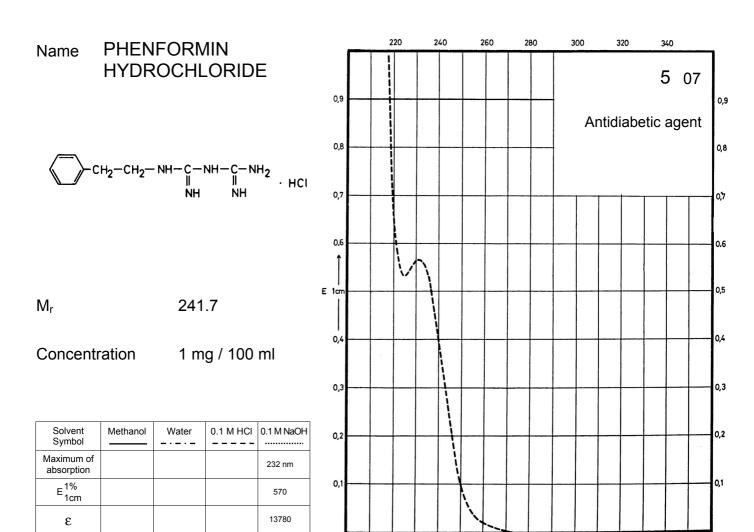


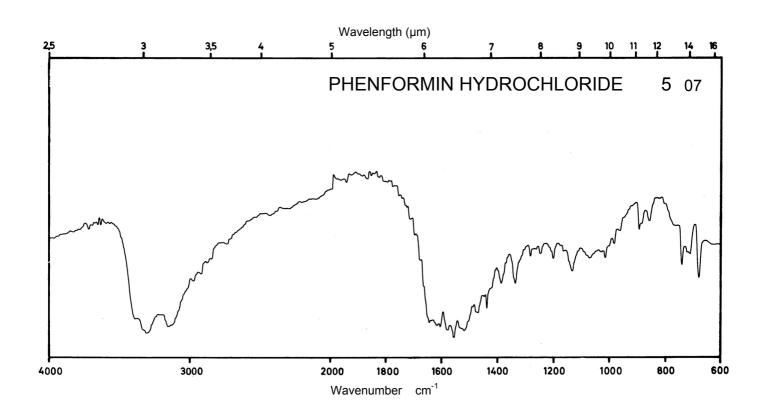


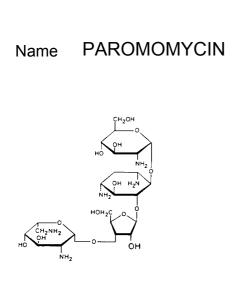



615.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

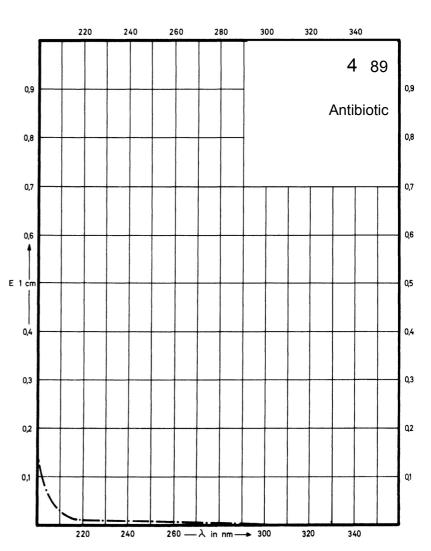


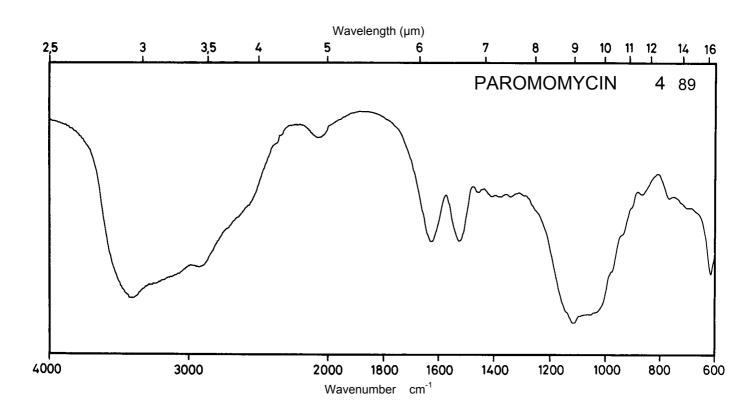

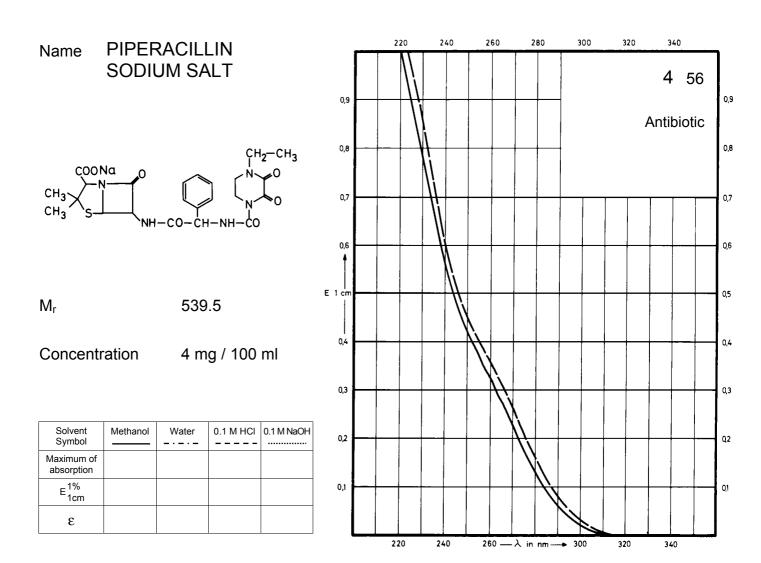

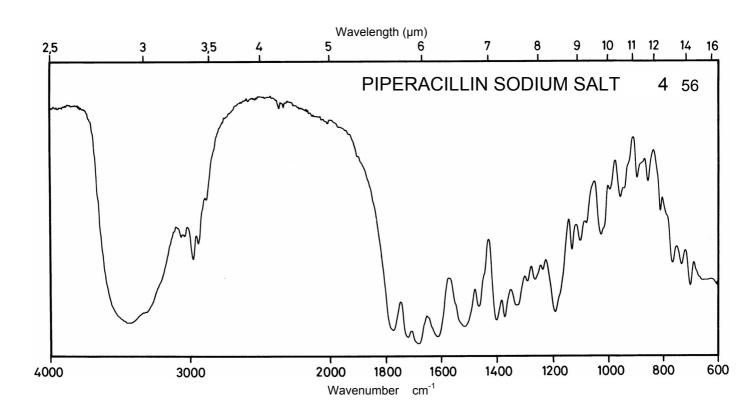




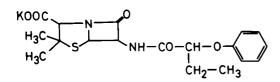



 M_{r}


416.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 271 nm	276 nm 269 nm	Decom- position observed	276 nm 269 nm
E ^{1%} 1cm	24.4 30.2	21.5 27.1		21.5 27.1
3	1020 1260	900 1130		900 1130

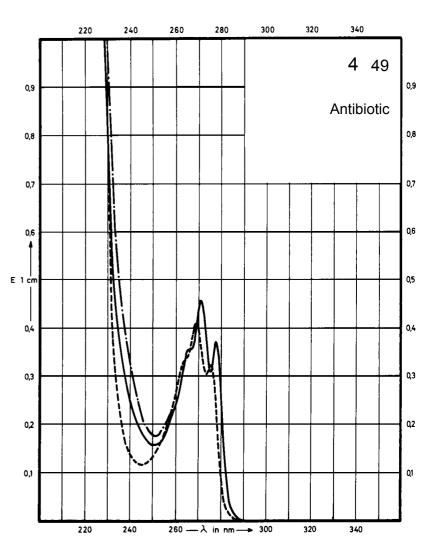


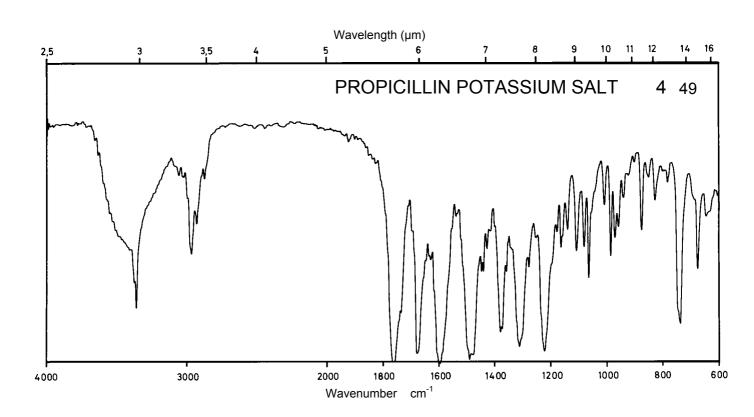

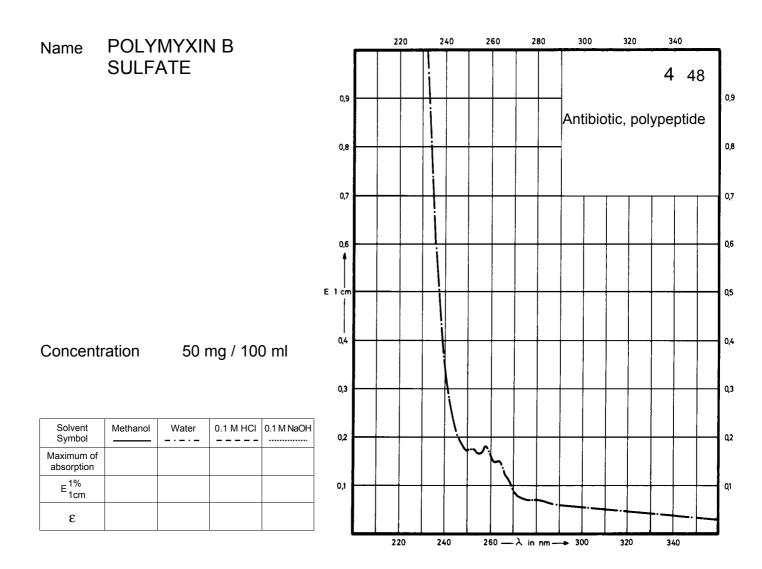

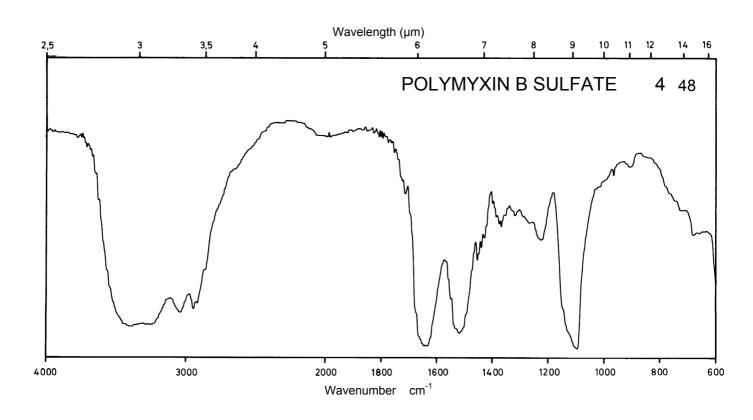




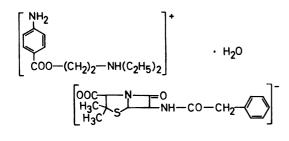



 M_r 588.7


Concentration 1 mg / 100 ml 13 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	295 nm		278 nm 226 nm	280 nm
E ^{1%} 1cm	356		51 274	248
ε	20960		3000 16130	14600

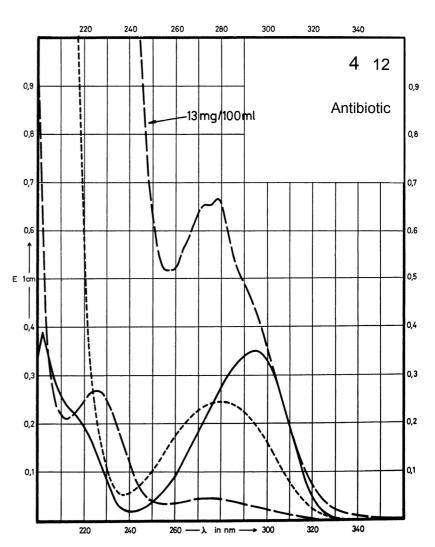

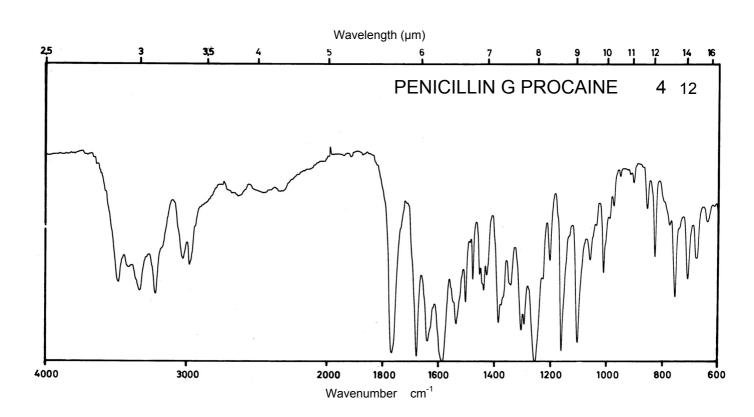


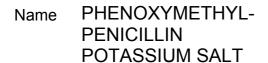


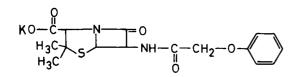



388.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 268 nm	275 nm 268 nm		273 nm 266 nm
E ^{1%} 1cm	29 35	25 31		26 32
ε	1130 1360	970 1200		1010 1240

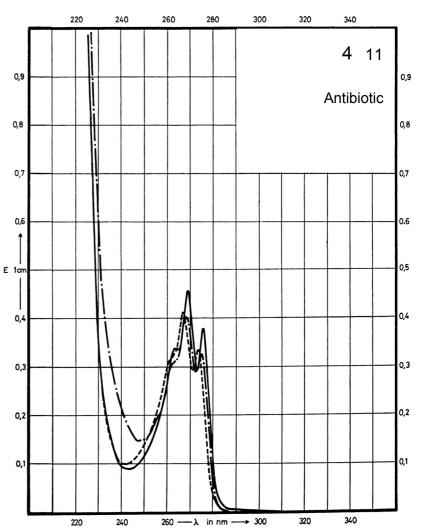


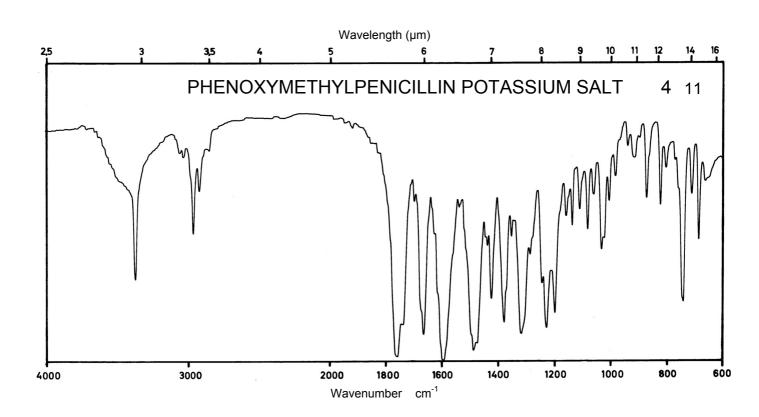

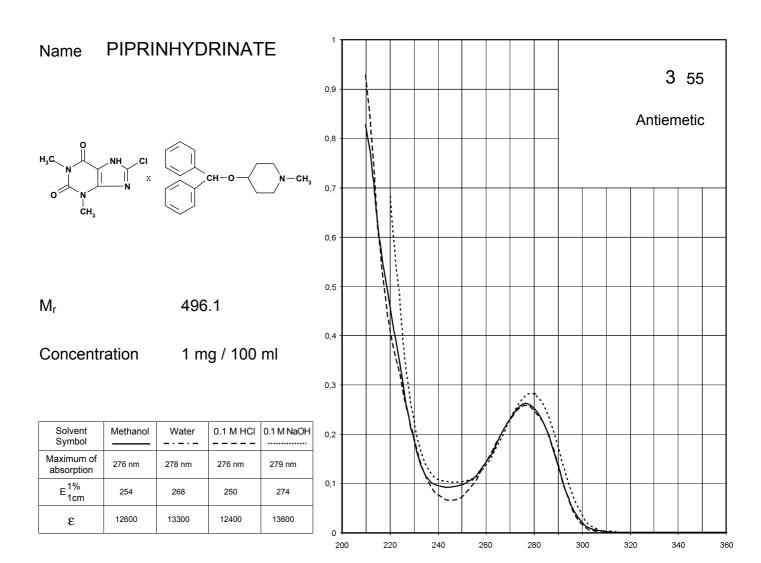

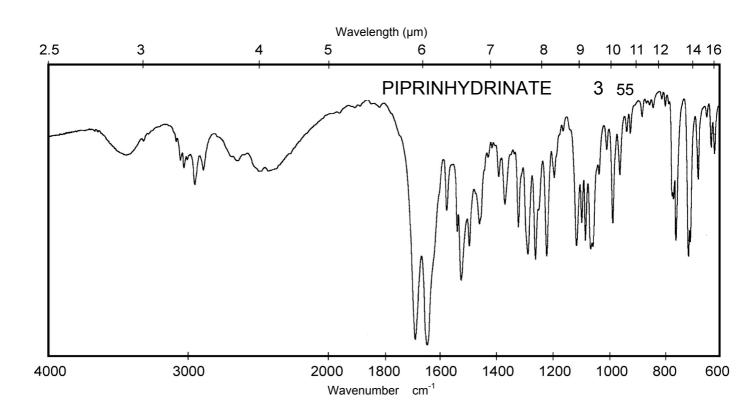

Name PIZOTIFEN MALEATE

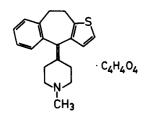


 M_{r}


411.5

Concentration 1.5


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	252 nm 231 nm		252 nm 229 nm	255 nm
E ^{1%} _{1cm}	249 412		245 407	225
3	10200 16900		10100 16800	9200

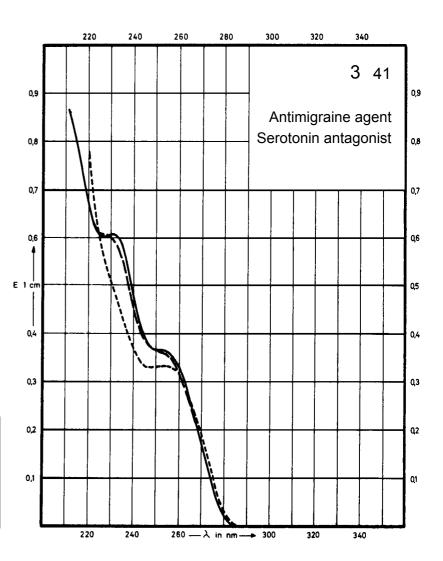

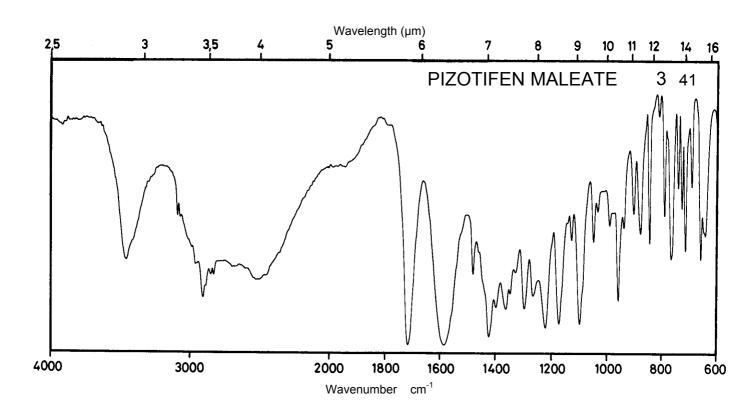





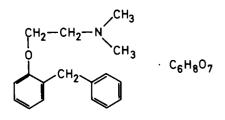

 M_r 447.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 269 nm		275 nm 269 nm	
E ^{1%} 1cm	37 43		38 45	
ε	1700 1900		1700 2000	

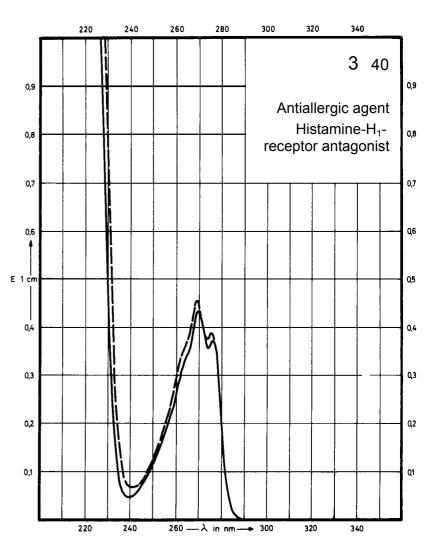

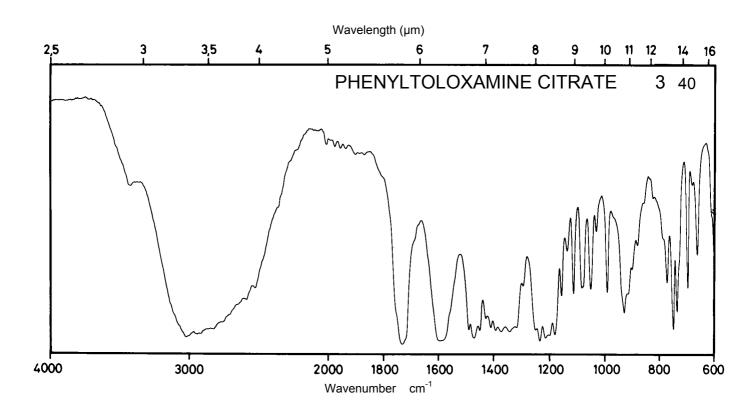

Name PROTHIPENDYL HYDROCHLORIDE


 M_{r}


321.9

Concentration

1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	325 nm 251 nm		372 nm 241 nm	320 nm 250 nm
E ^{1%} 1cm	143 920		159 740	137 760
ε	4600 29610		5120 23820	4410 24460

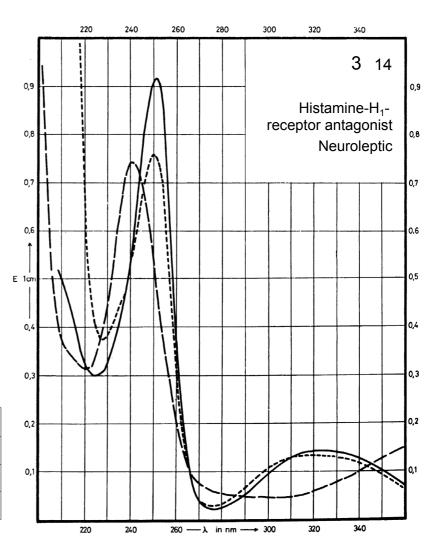

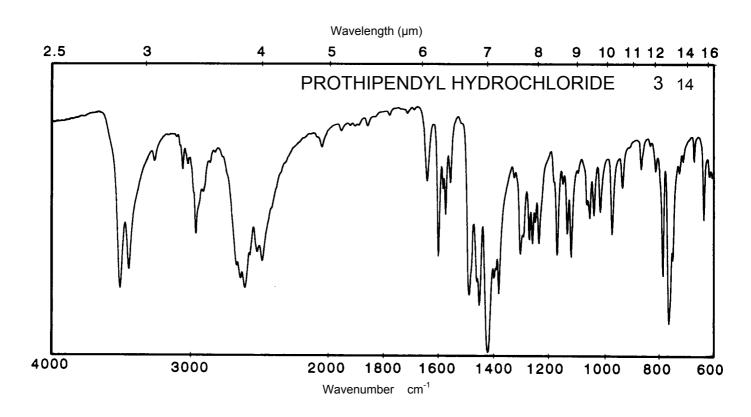



Name PHENIRAMINE 4-AMINOSALICYLATE

 M_r


393.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	301 nm 264 nm		299 nm 266 nm	298 nm 262 nm
E ^{1%} 1cm	268 409		135 225	217 446
3	10550 16090		5310 8850	8540 17550

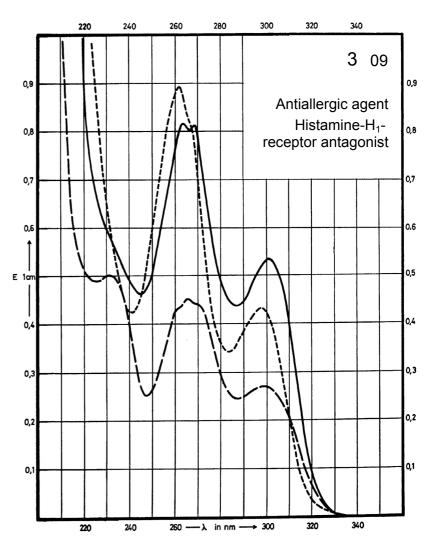

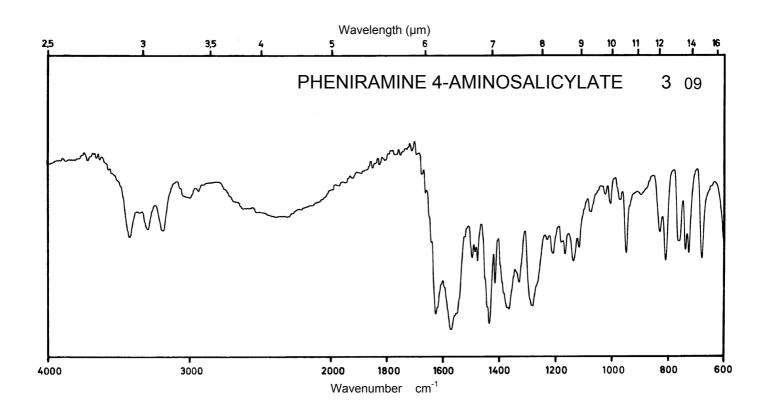





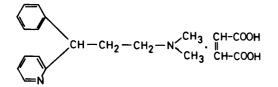


356.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm		264 nm	263 nm
E ^{1%} 1cm	138		210	138
3	4920		7480	4920

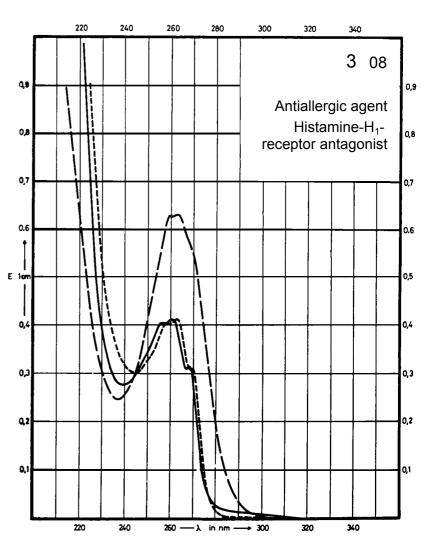

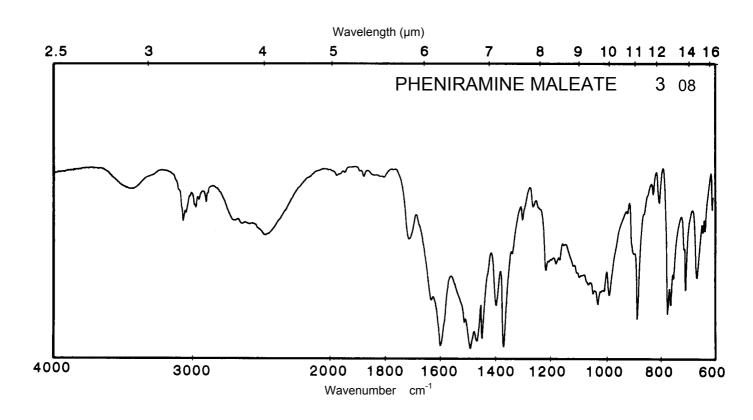





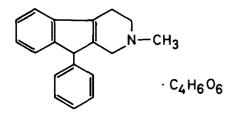


 M_{r} 411.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	260 nm		260 nm	262 nm
E ^{1%} 1cm	233		229	248
3	9590		9420	10210

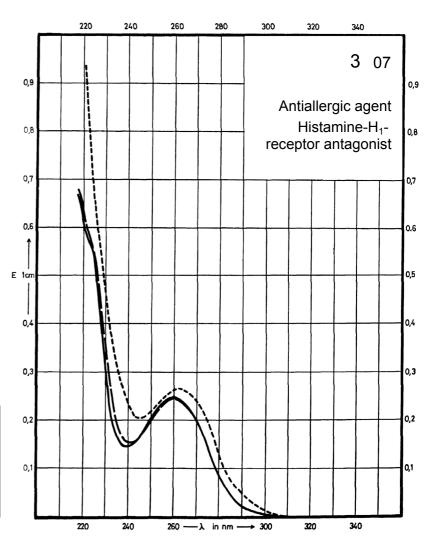

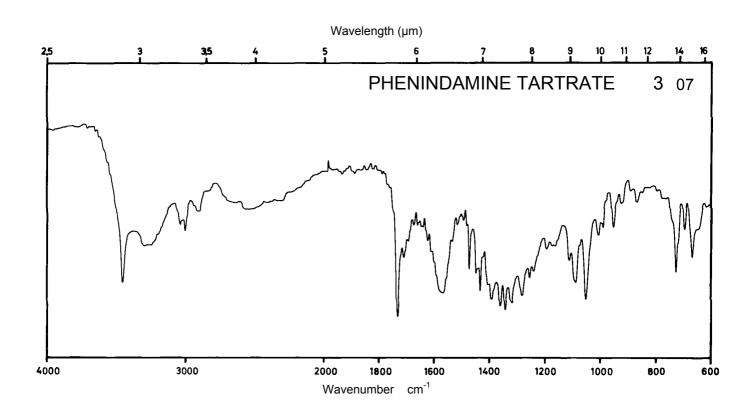





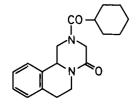


312.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 264 nm			
E ^{1%} 1cm	10.6 11.4			
ε	330 355			

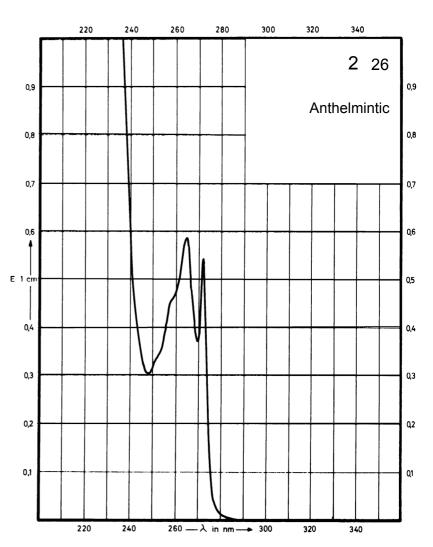

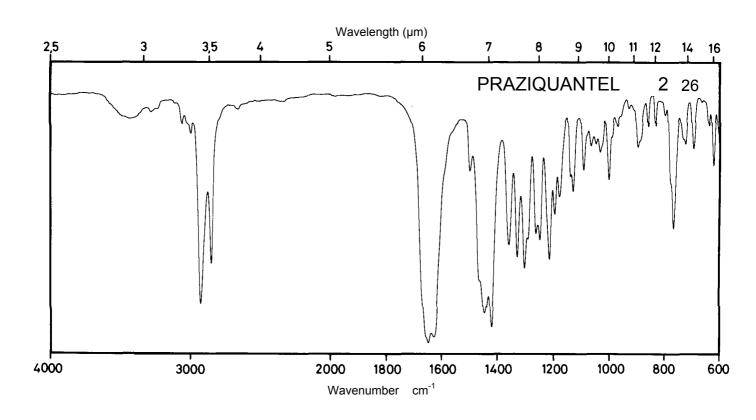





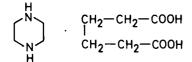


232.3


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

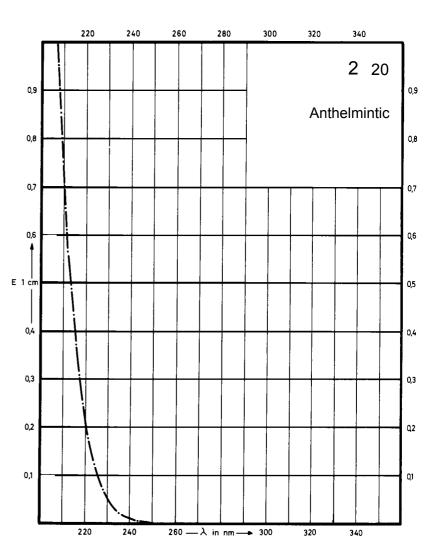


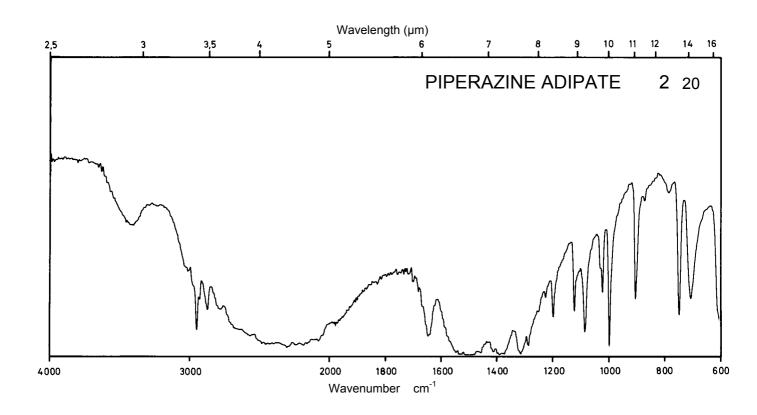

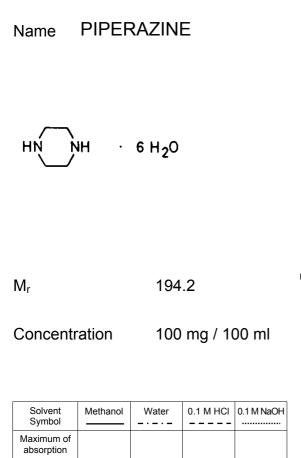


ε

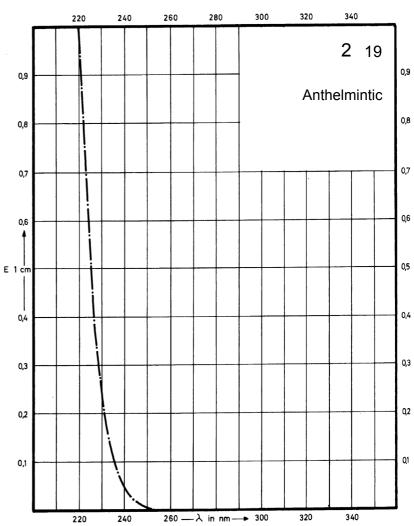

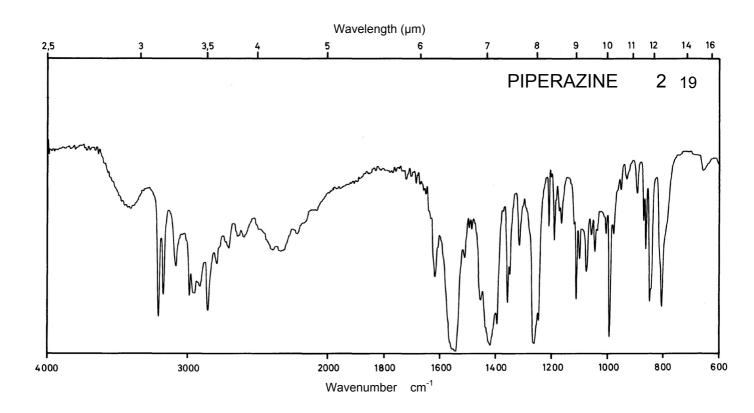



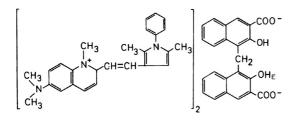


1151.4


Concentration


0.5 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI*	0.1 M NaOH*
Maximum of absorption	505 nm 356 nm 235 nm		500 nm 356 nm 236 nm	500 nm 356 nm 235 nm
E ^{1%} 1cm	750 390 1315		695 355 1350	695 375 1320
3	86400 44900 151400		80000 40900 155400	80000 43200 152000

^{* 1} M HCl + methanol (1+9)

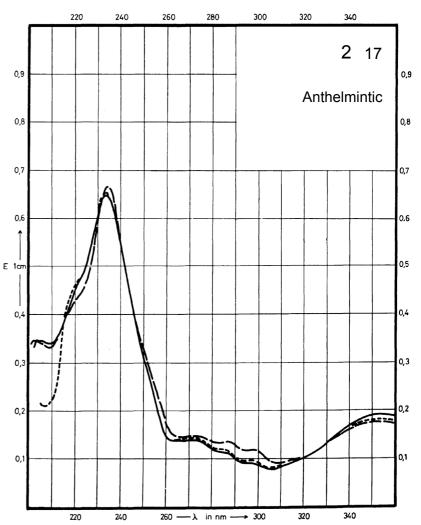

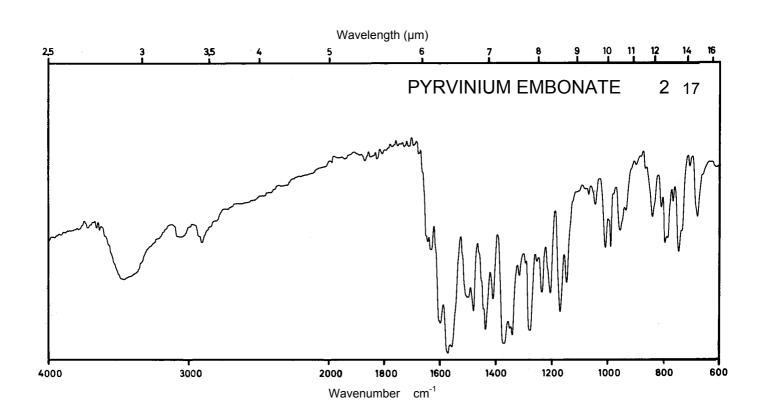

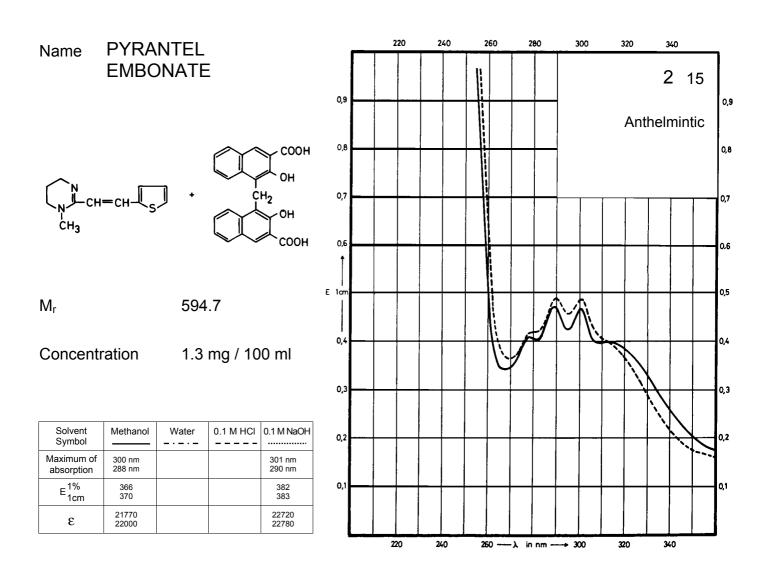

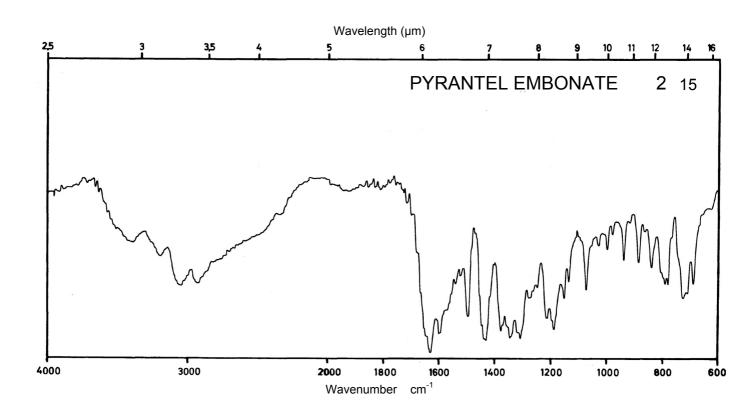



Name PARBENDAZOLE

 M_{r}


247.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	292 nm 285 nm 246 nm		287 nm 281 nm	303 nm 297 nm 254 nm
E ^{1%} 1cm	651 579 493		737 639	815 816 462
3	16100 14320 12190		18230 15800	20150 20180 11430

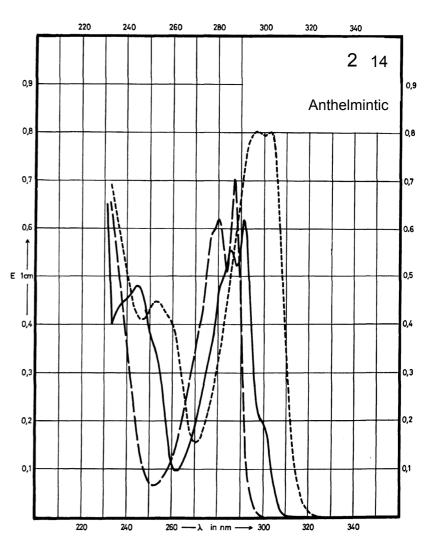

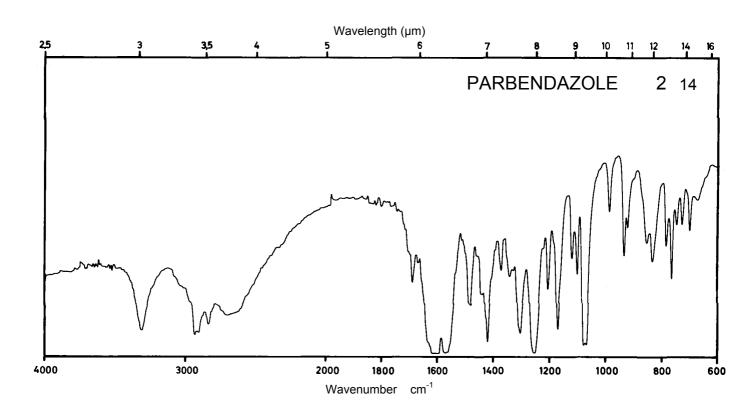



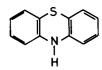



199.3

Concentration


0.25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 252 nm		250 nm	249 nm
E ^{1%} 1cm	237 2050		1350	1460
ε	4720 40860		26910	29100

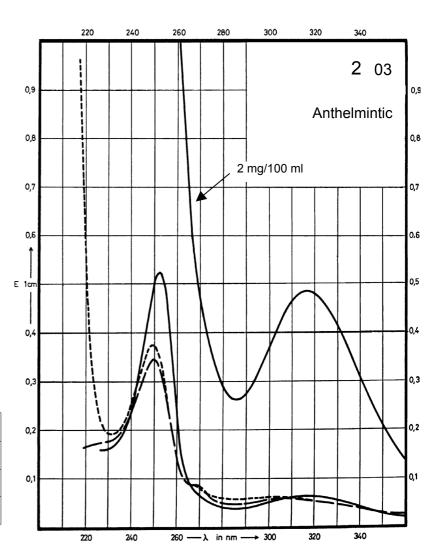

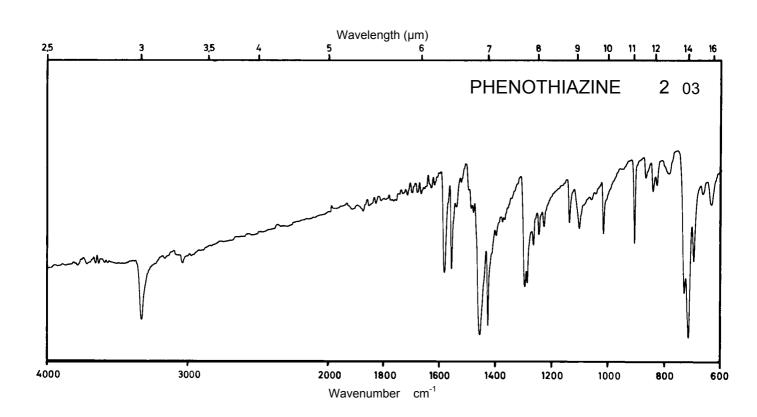




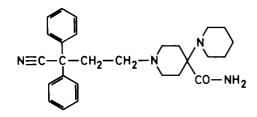


 M_r 430.6


Solvent Symbol	Methanol	Water	0.1 M HCI*	0.1 M NaOH
Maximum of absorption			264 nm 258 nm 253 nm	
E ^{1%} 1cm			8.0 10.2 8.7	
ε			344 438 375	

^{* 1} M HCl + Propanol-2 (1+9)

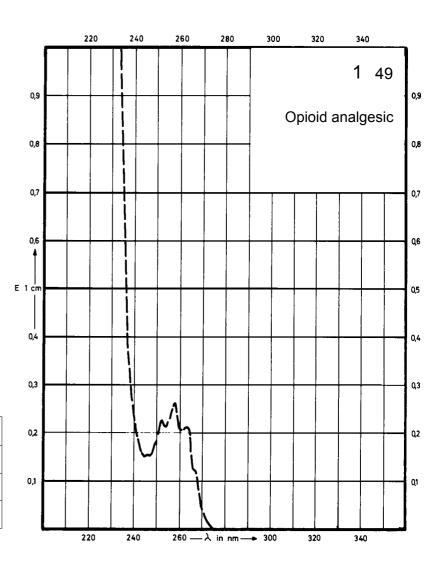

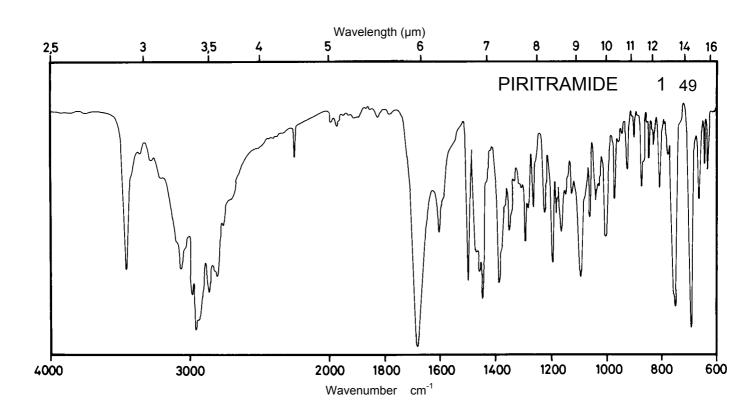





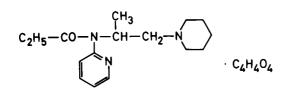


391.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm		258 nm	258 nm
E ^{1%} 1cm	112		107	119
ε	4400		4200	4700

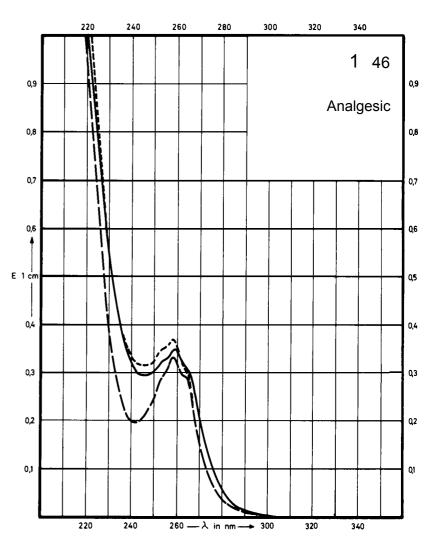

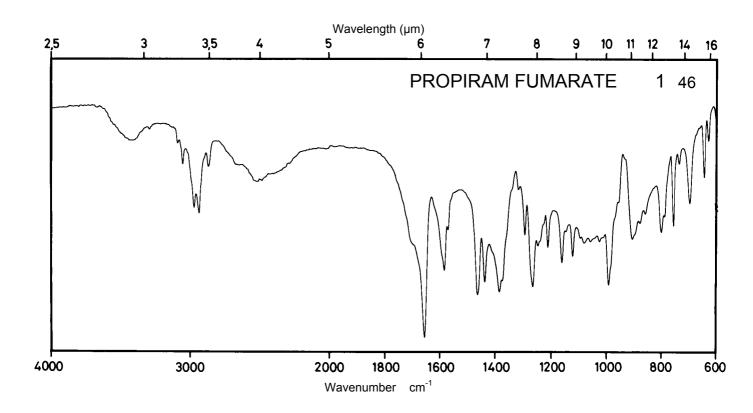




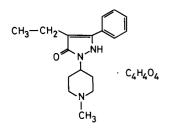


 M_{r}


401.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	247 nm		249 nm	247 nm
E ^{1%} 1cm	290		322	333
3	11600		12900	13400

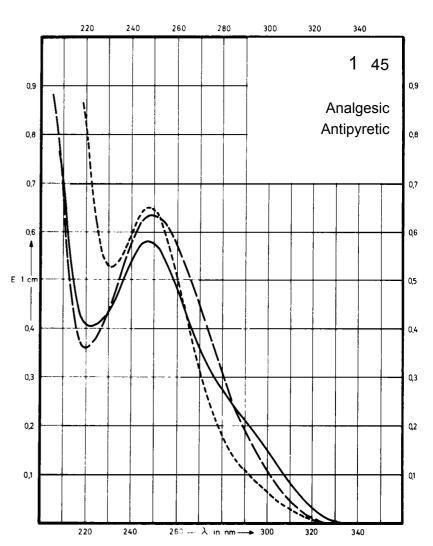

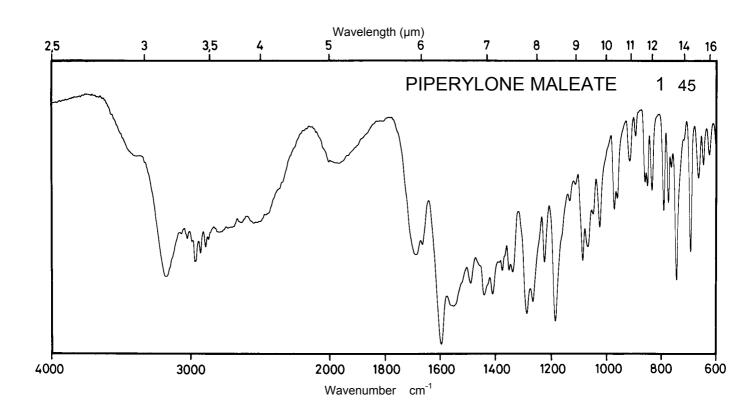





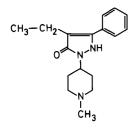


 M_r 285.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	247 nm		249 nm	247 nm
E ^{1%} 1cm	457		468	485
ε	13000		13400	13800

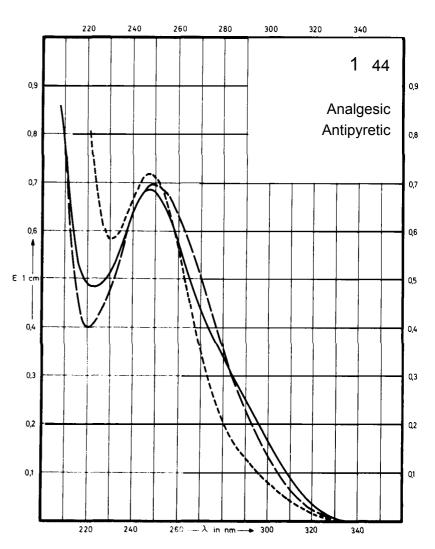

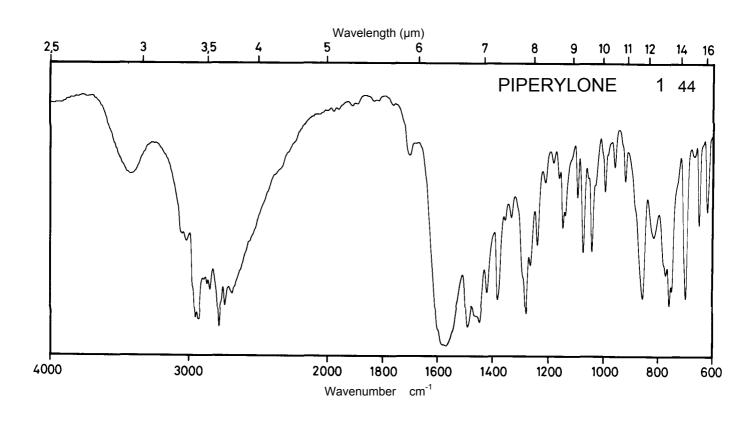


Name PHENAZONE COMP. WITH SALICYLIC ACID

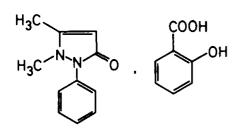


 M_{r} 326.3


Concentration 0.9 mg / 100 ml 1.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm 238 nm		303 nm 234 nm	296 nm
E ^{1%} 1cm	294 489		111 585	110
3	9600 15950		3630 19100	3580

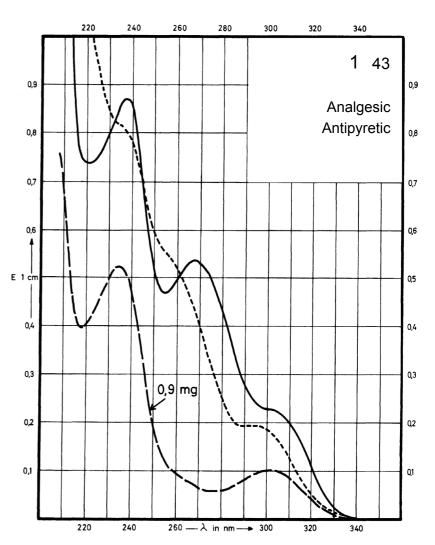

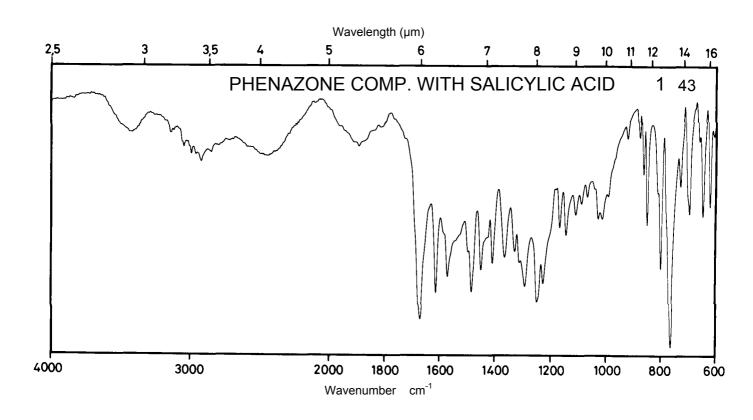





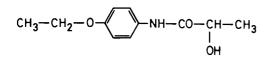


 M_r 209.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	249 nm		244 nm	245 nm
E ^{1%} 1cm	740		535	557
3	15500		11200	11700

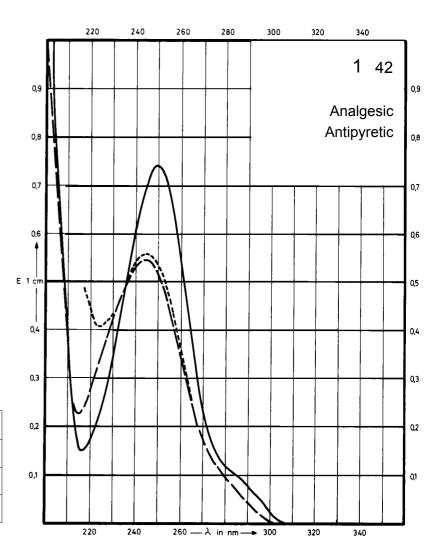

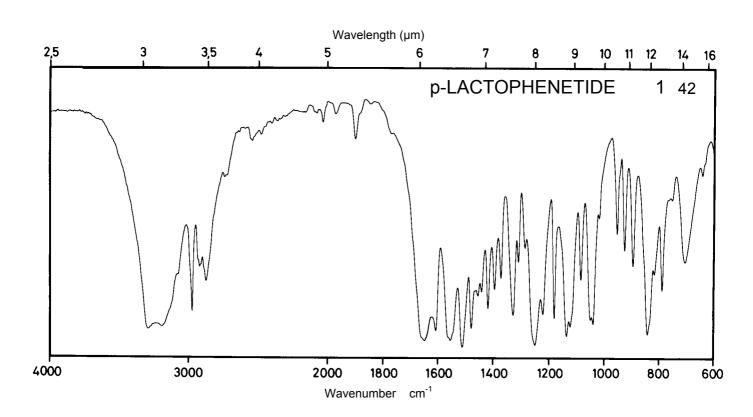




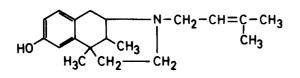


285.4


Concentration 2 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		278 nm	298 nm 239 nm
E ^{1%} 1cm	82		70	107 326
3	2330		1990	3070 9310

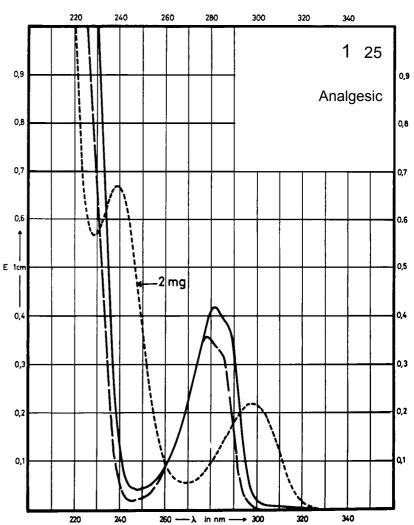

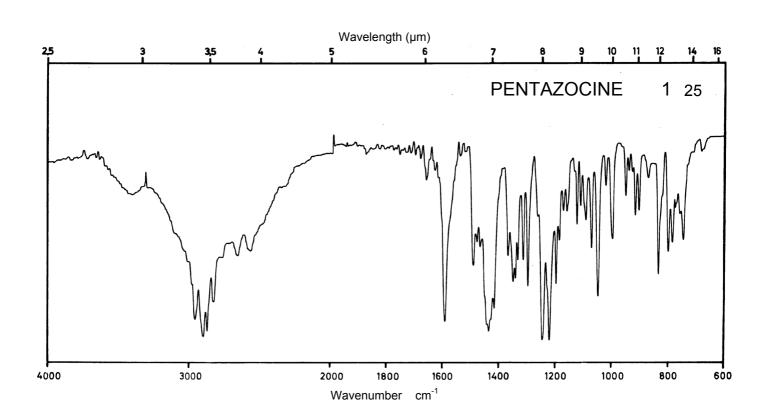




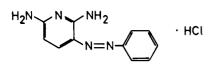


249.7


Concentration

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	395 nm 238 nm		425 nm 239 nm	430 nm 276 nm
E ^{1%} 1cm	935 415		755 397	905 435
ε	23350 10360		18850 9910	22600 10860

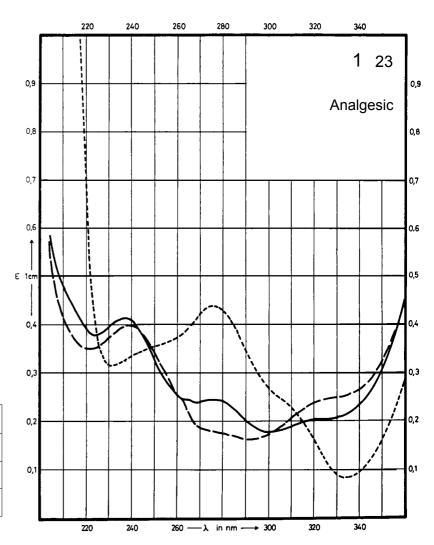

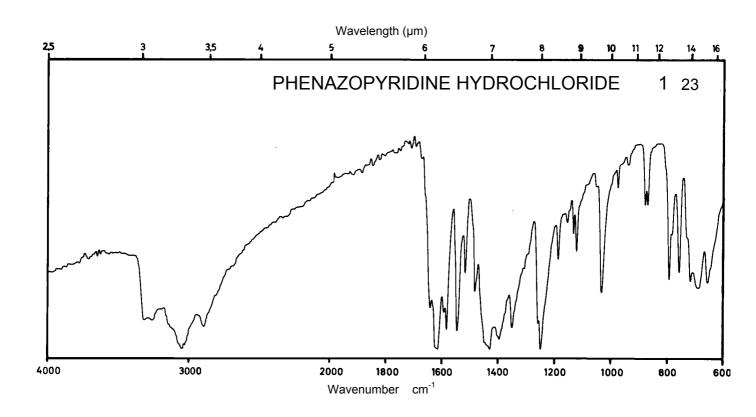




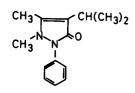


 M_r


230.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	275 nm 246 nm		240 nm	265 nm 245 nm
E ^{1%} 1cm	420 425		400	390 385
ε	9670 9790		9210	8980 8870

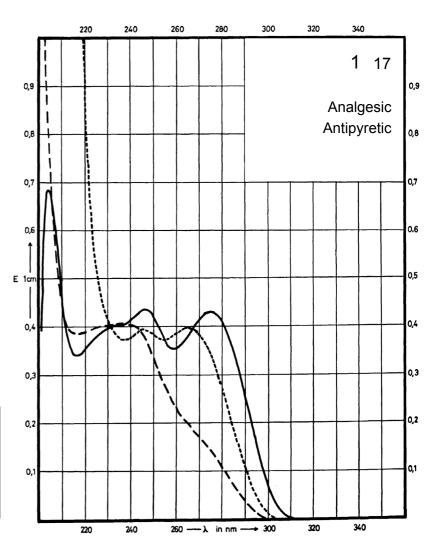

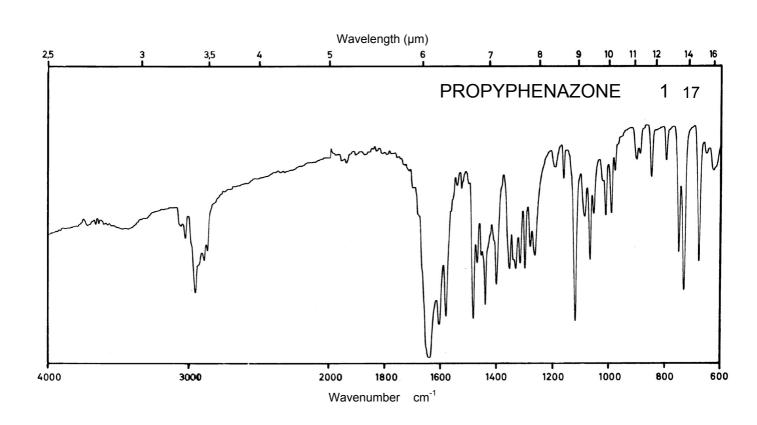




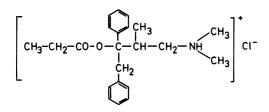


 M_{r}

375.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm		257 nm	
E ^{1%} 1cm	10.5		10.9	
ε	390		410	

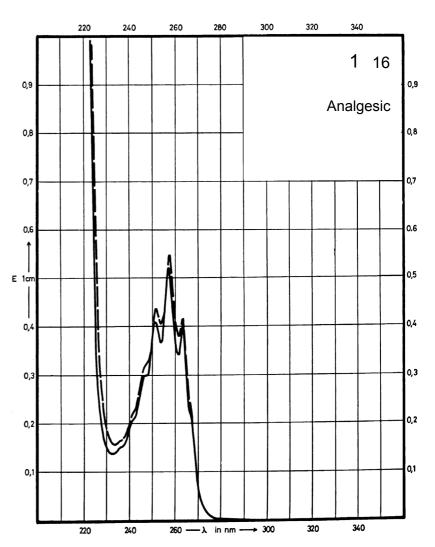

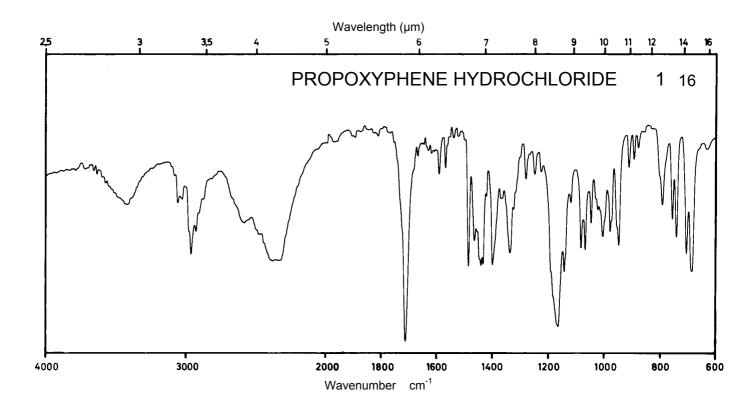




H₃C N N O


 M_{r}

188.2

Concentration

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	267 nm 243 nm	256 nm 243 nm	231 nm	256 nm 243 nm
E ^{1%} 1cm	492 504	471 492	597	471 494
ε	9260 9480	8860 9260	11230	8860 9260

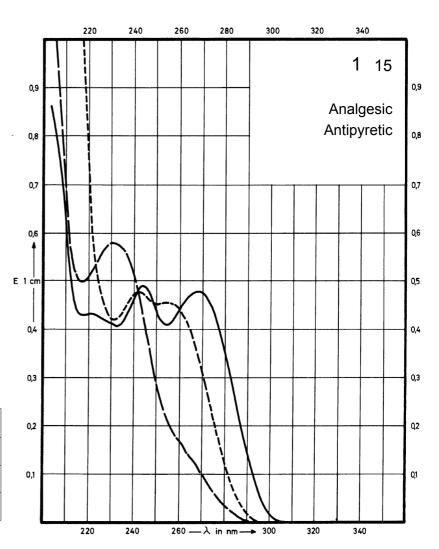

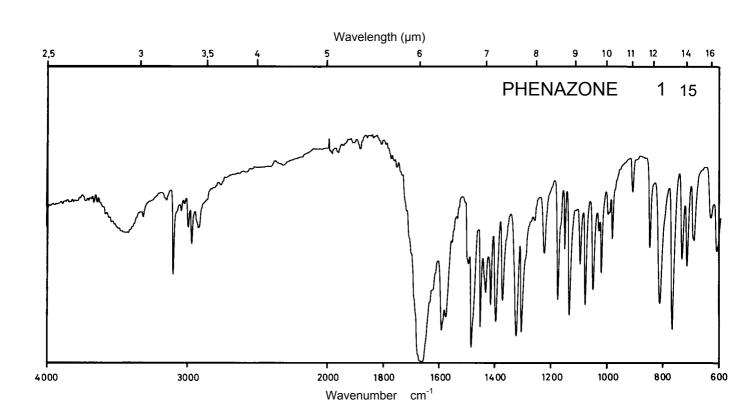

Name PHENACETIN


NH-CO-CH₃


 M_{r}

179.2


Concentration 1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
aximum of bsorption	249 nm		244 nm	245 nm
E ^{1%} _{1cm}	864		647	657
ε	15480		11590	11770

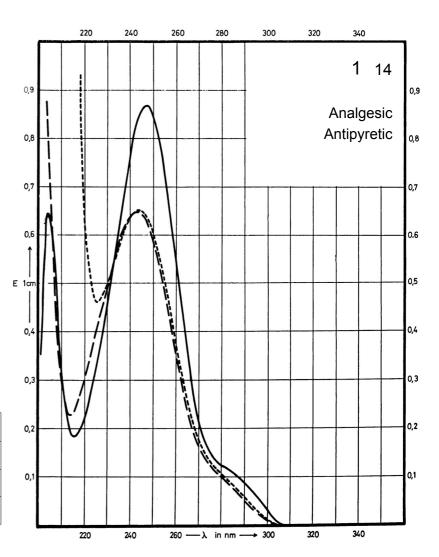

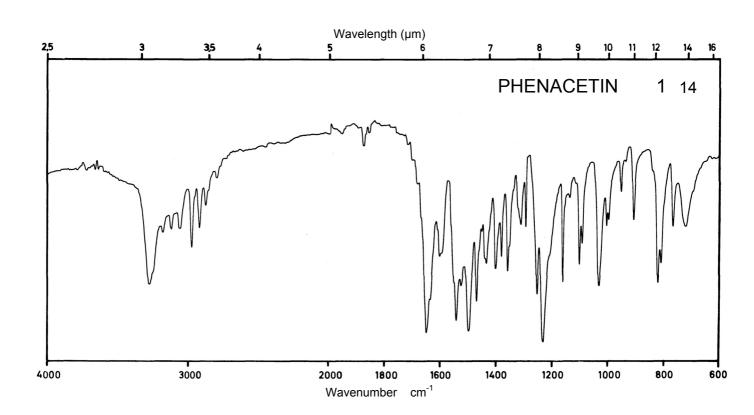




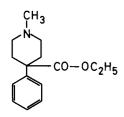


247.3


Concentration

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 259 nm 253 nm			
E ^{1%} 1cm	6.4 8.2 7.0			
ε	157 204 172			

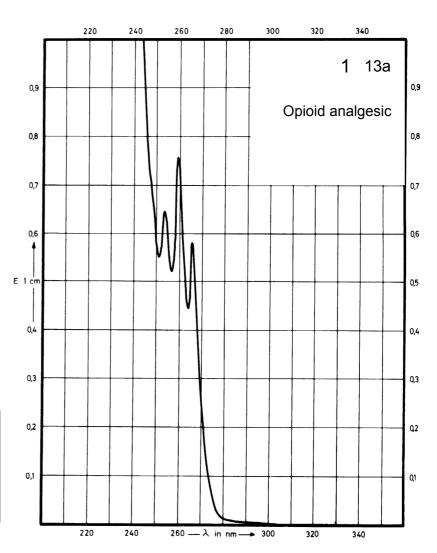

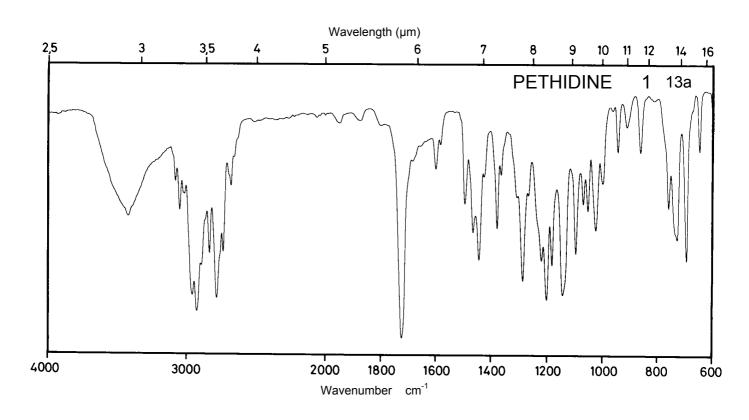


Name PETHIDINE HYDROCHLORIDE

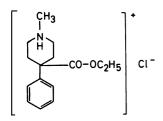


 M_{r}

283.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 257 nm 251 nm		263 nm 257 nm 251 nm	263 nm 257 nm 251 nm
E ^{1%} 1cm	5.5 7.2 6.0		6.0 7.5 6.2	5.7 7.4 6.3
ε	156 204 170		170 213 176	162 210 179

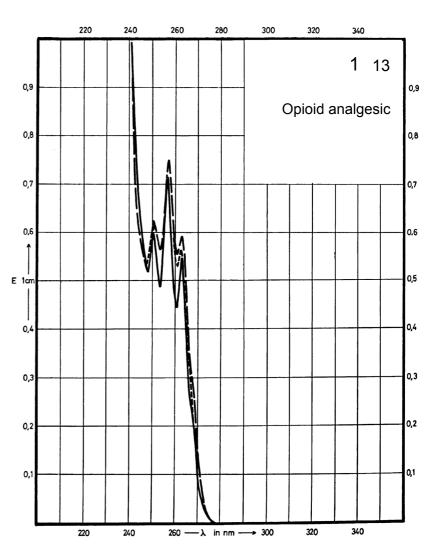

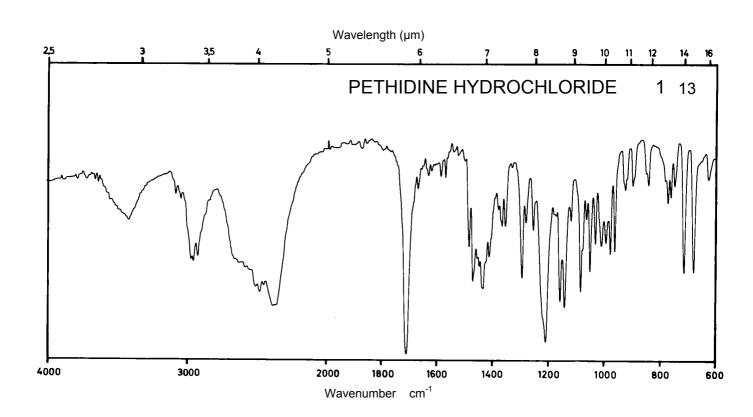



Name PARACETAMOL

OH NH-CO-CH₃


 M_r 151.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	247 nm		240 nm	255 nm
E ^{1%} 1cm	850		642	710
3	12850		9710	10740

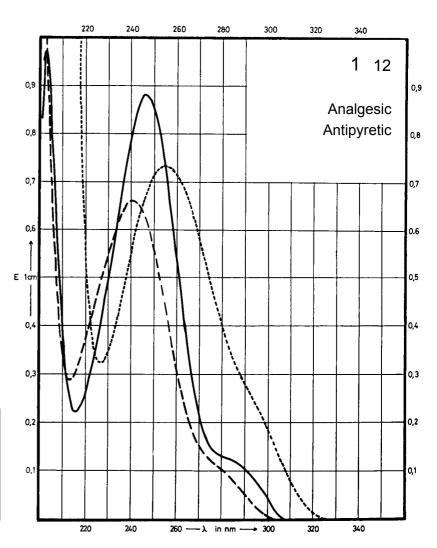

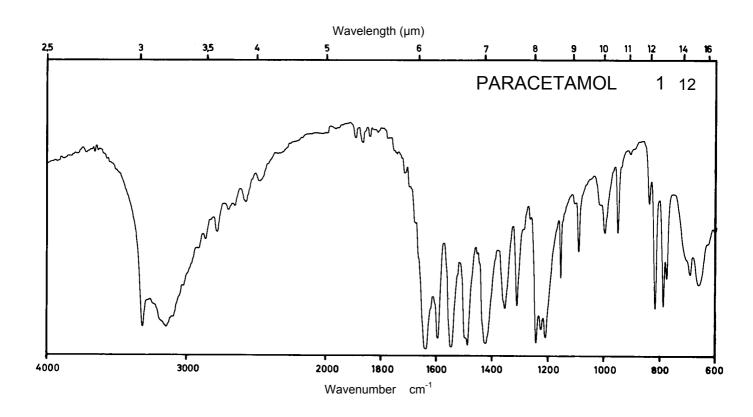





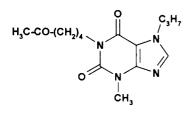


306.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm	274 nm	273 nm	274 nm
E ^{1%} 1cm	308	322	313	320
ε	9400	9900	9600	9800

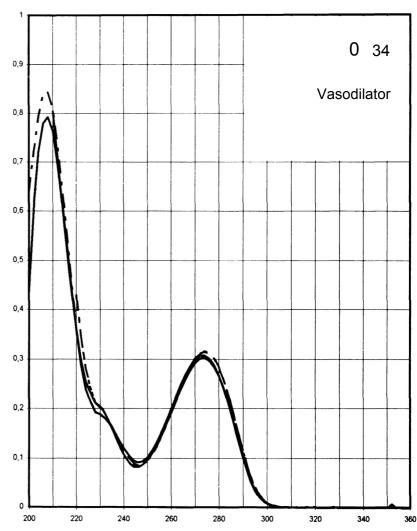

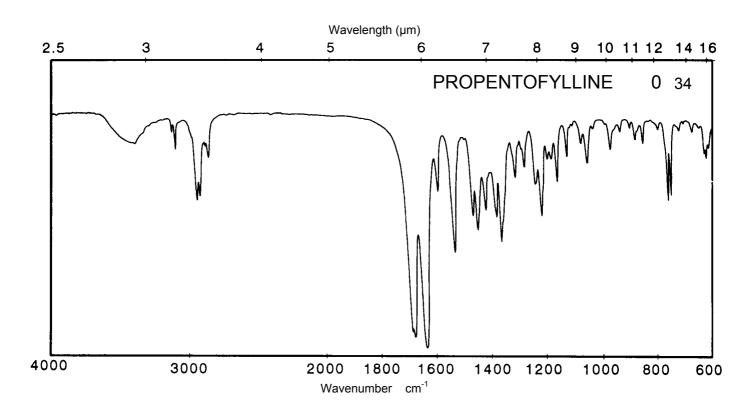





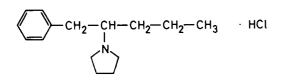


253.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm		263 nm 257 nm 252 nm	
E ^{1%} 1cm	5.61 7.20 5.57		5.70 7.37 5.96	
ε	140 180 140		145 190 150	

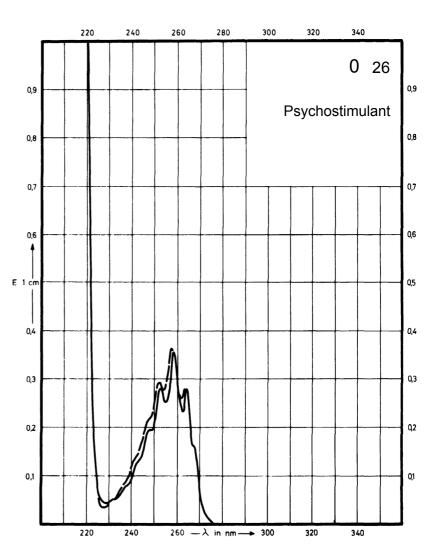

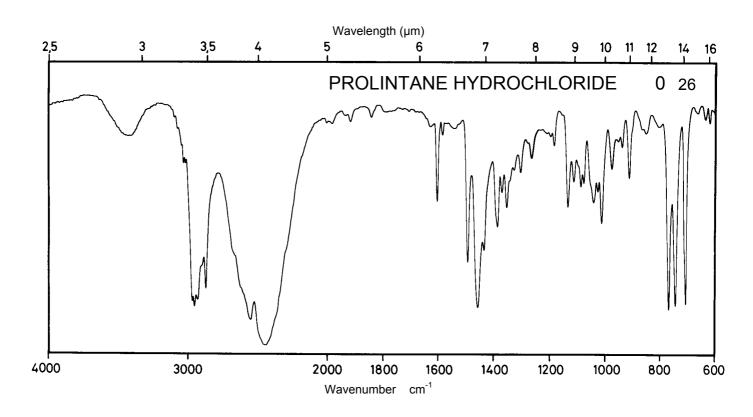





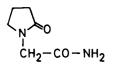


 M_r 142.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
ε				

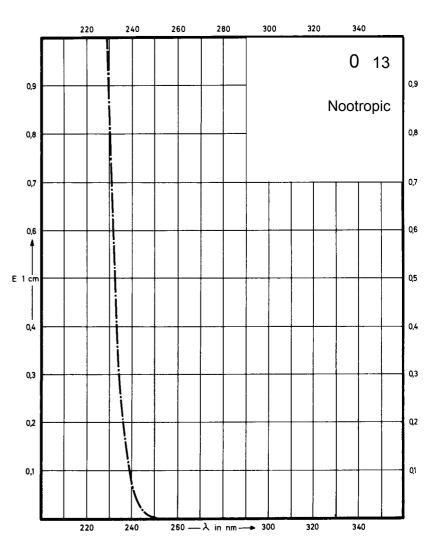

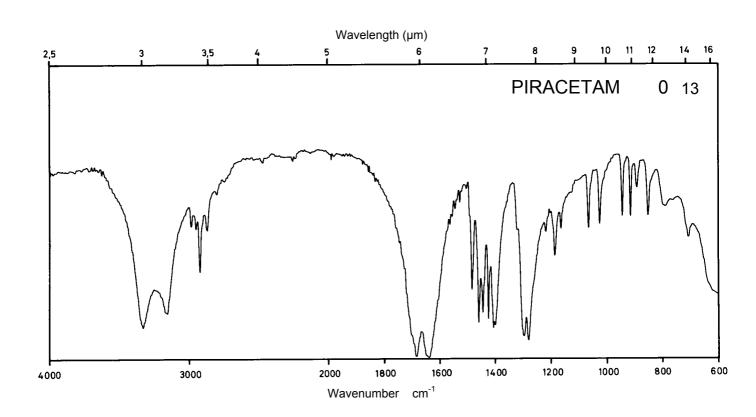



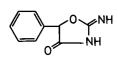

 M_r

Concentration 10 mg / 100 ml 40 mg / 100 ml


176.2


 Solvent Symbol
 Methanol --- Water --- 0.1 M HCI --- 0.1 M NaOH 0.1 M NaOH

 Maximum of absorption
 268 nm 262 nm 262 nm 256 nm
 268 nm 261 nm


 E 1% 1cm
 11.1 19.2
 12.5 19.9 22.6
 7.9 12.9

 E
 195 340
 350 350 398
 140 227

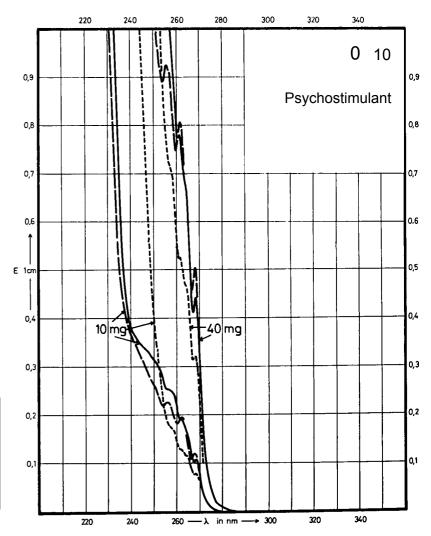

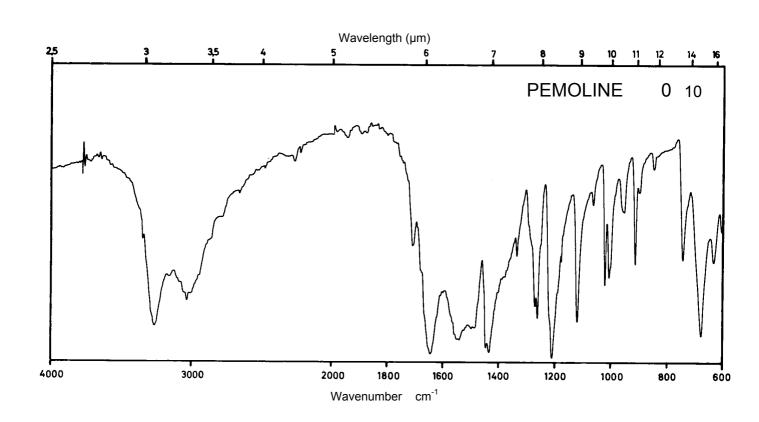





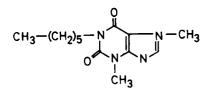


264.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		273 nm	273 nm
E ^{1%} 1cm	367		373	382
ε	9700		9860	10100

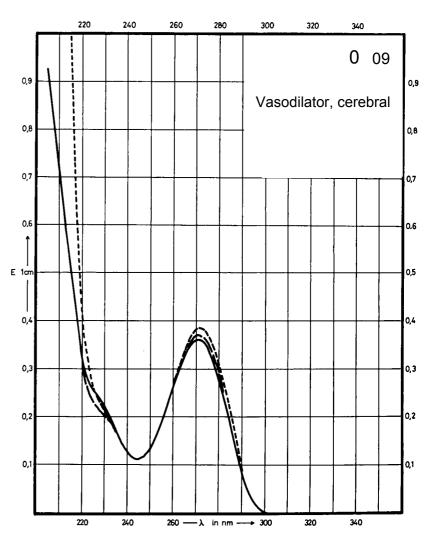

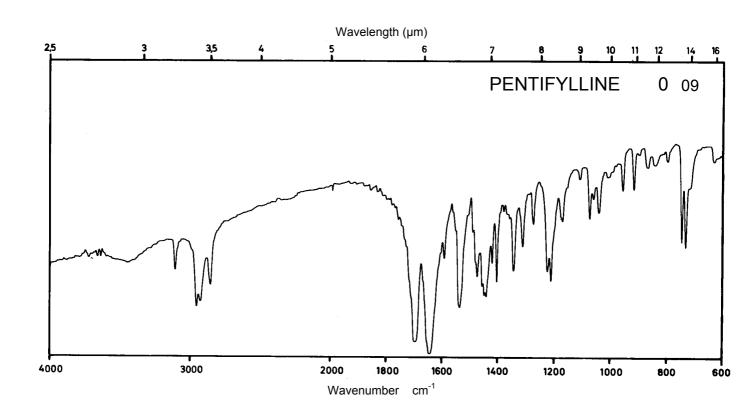





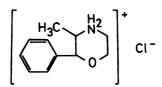


213.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	255 nm			
E ^{1%} 1cm	8.6			
3	180			

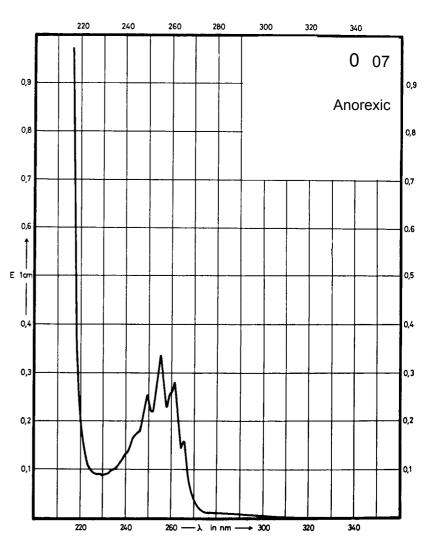

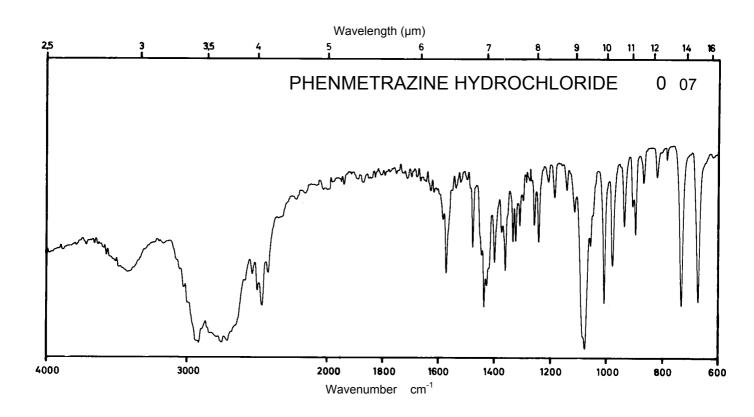




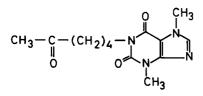


 M_r

278.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		274 nm	272 nm
E ^{1%} 1cm	342		359	365
ε	9520		9990	10160

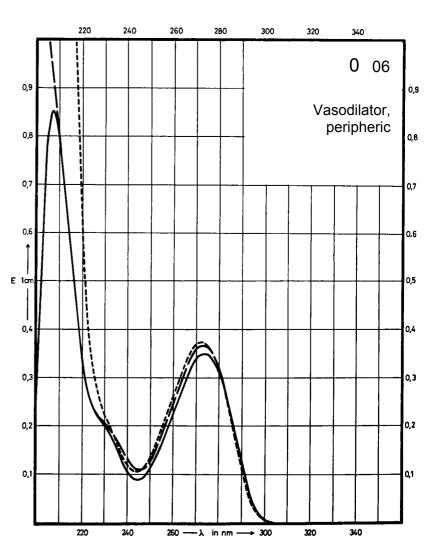

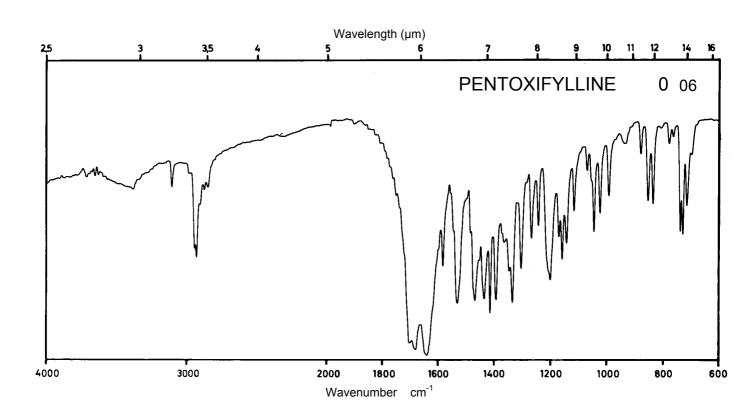







138.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

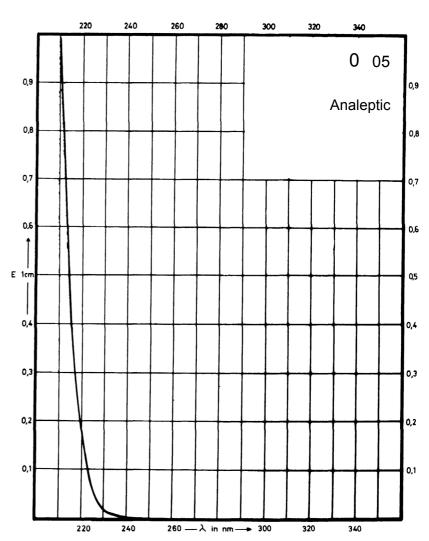

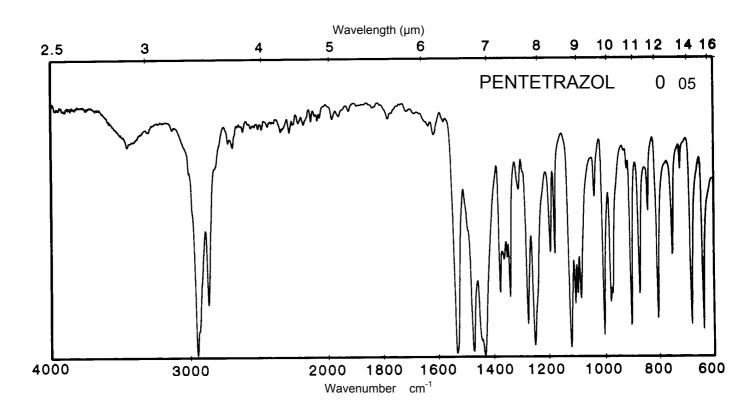



Name PHENIRAMINE

 M_{r}

240.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm 269 nm		265 nm	262 nm 269 nm
E ^{1%} 1cm	182 143		302	191 145
3	4380 3440		7260	4590 3490

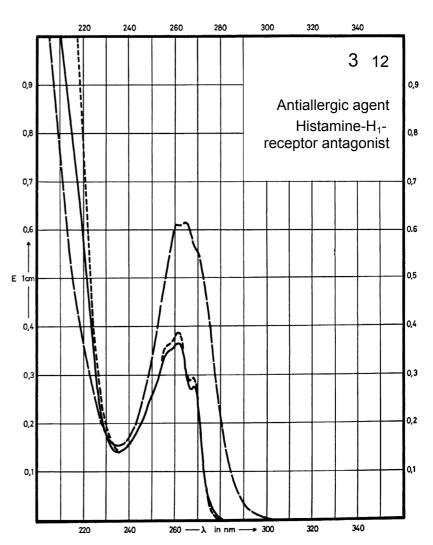

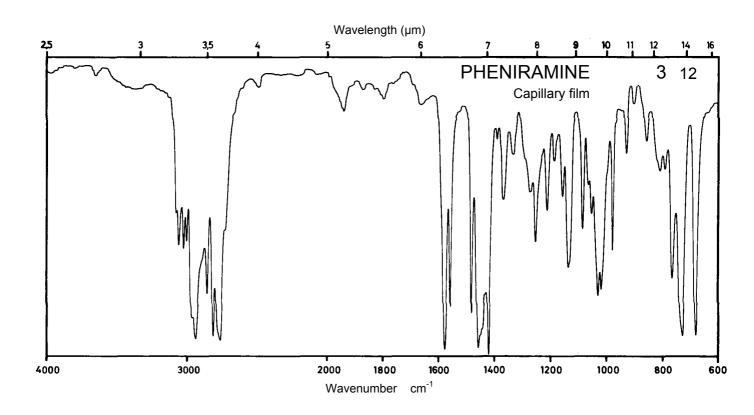







114.2

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm	252 nm	252 nm	244 nm
E ^{1%} 1cm	1360	1350	1360	1160
ε	15500	15400	15500	13200

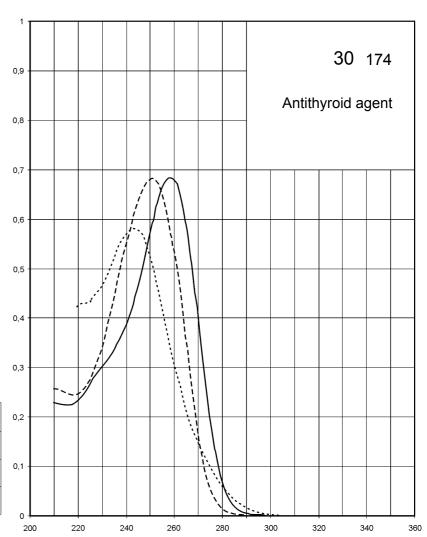

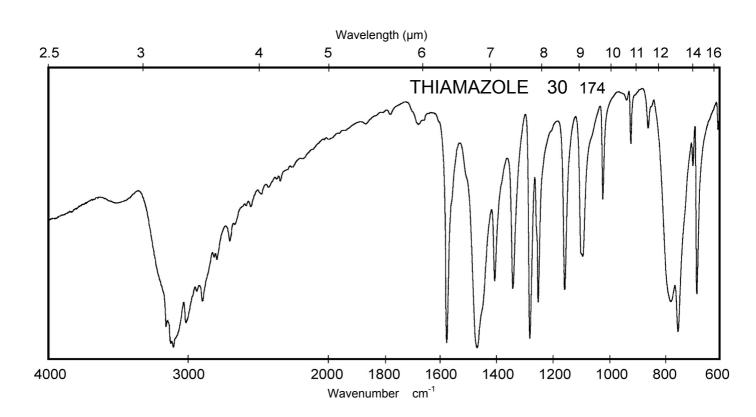





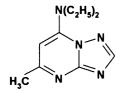


 M_r 205.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	306 nm 222 nm	Decom- position observed	301 nm 221 nm	307 nm 270 nm
E ^{1%} 1cm	856 906		908 888	881 254
3	17500 18600		18600 18200	18100 5200

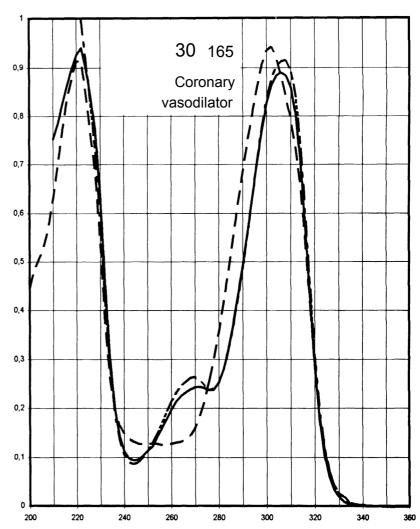

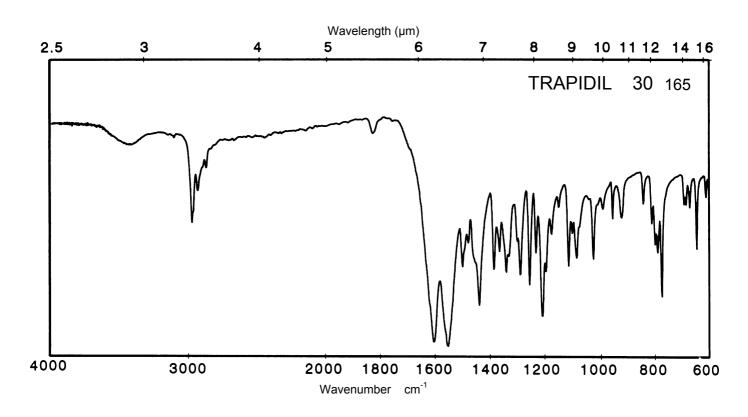




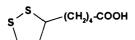



206.3

Concentration 1 mg / 100 ml 50 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	334 nm			
E ^{1%} 1cm	6.9			
ε	140			

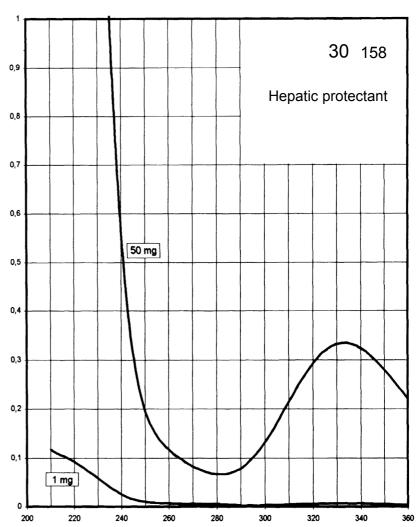


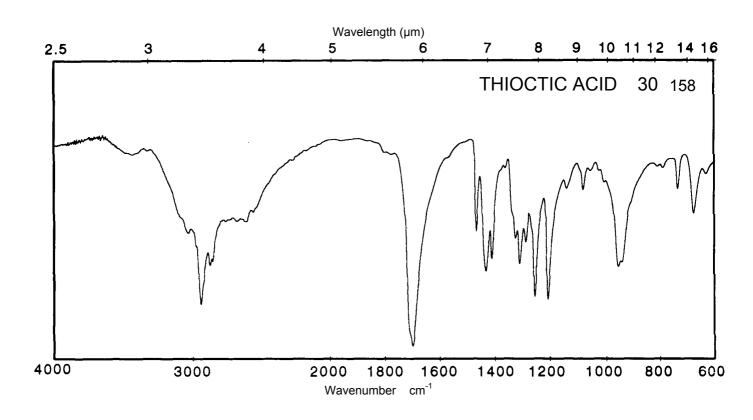

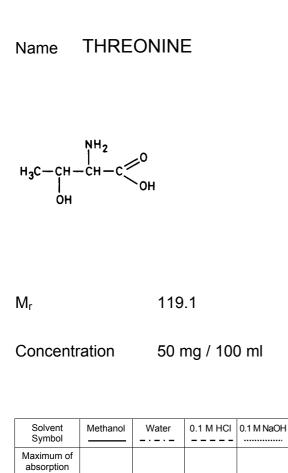


ε

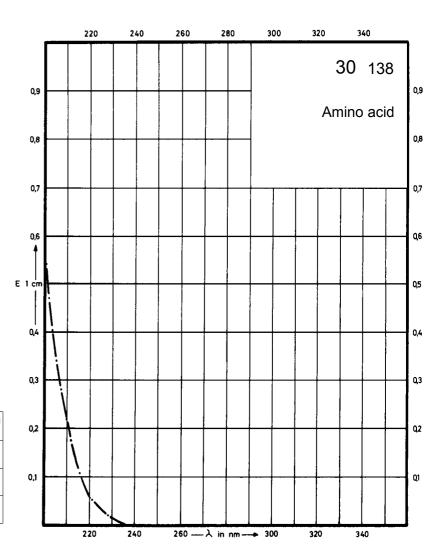

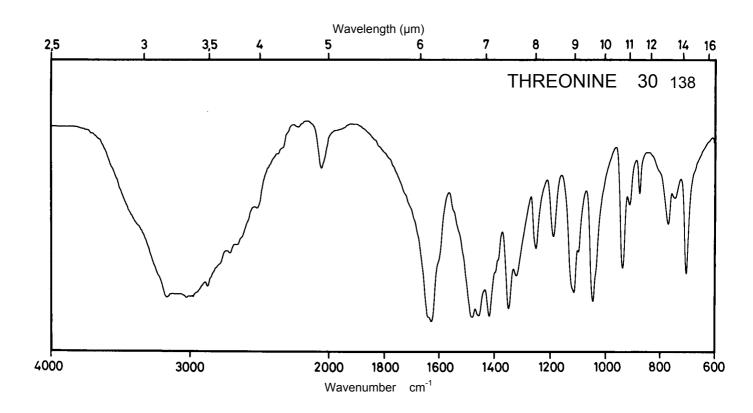





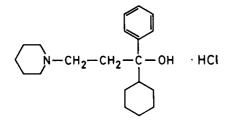


337.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 258 nm 252 nm			
E ^{1%} 1cm	4.33 5.92 5.00			
ε	146 200 170			

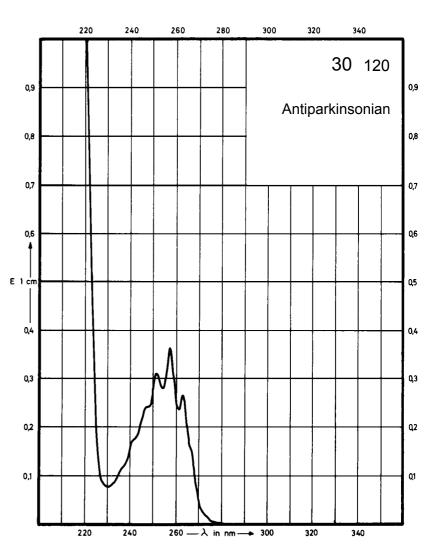

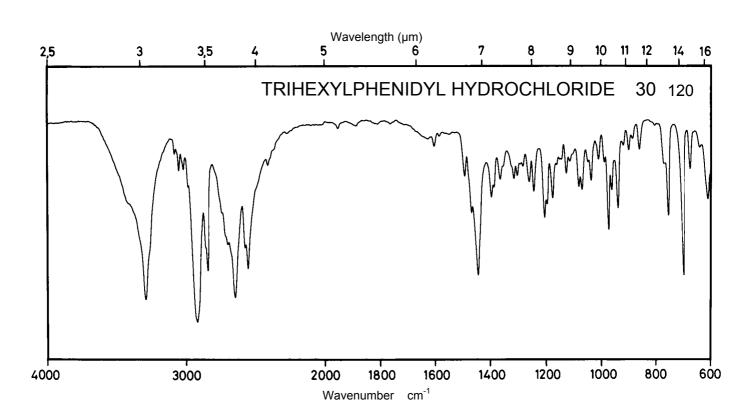





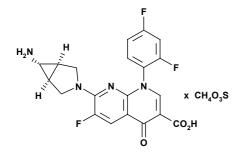



512.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	342 nm	344 nm	348 nm	346 nm
	270 nm	267 nm	272 nm	264 nm
E ^{1%}	397	386	365	418
1cm	877	736	837	683
3	20300	19800	18700	21400
	44900	37700	42900	35000

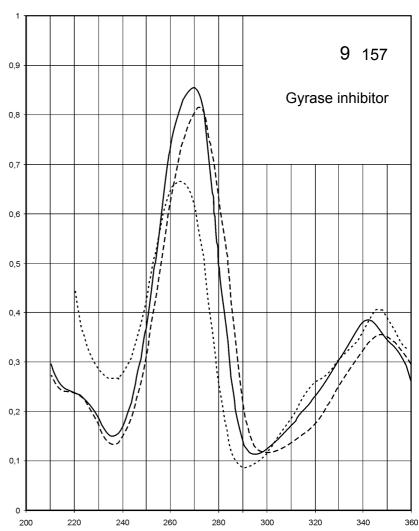


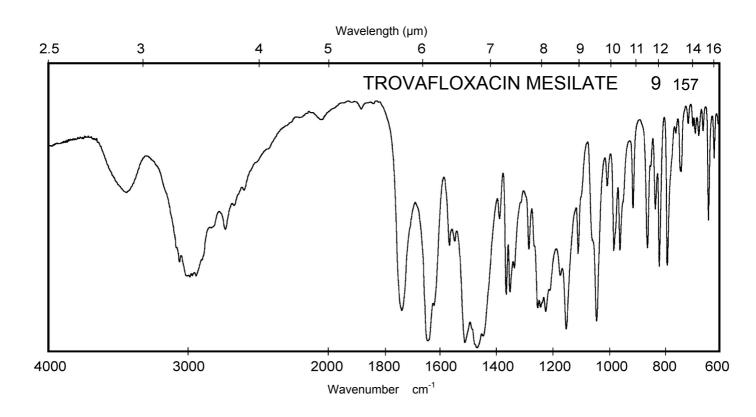

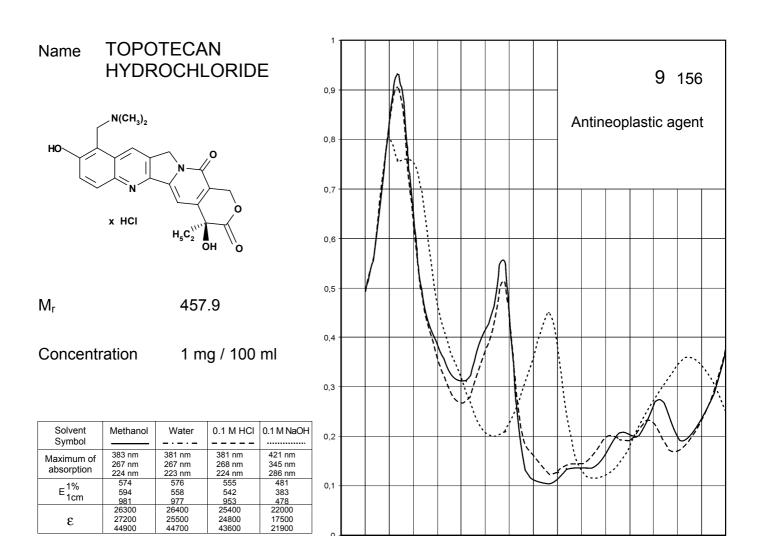

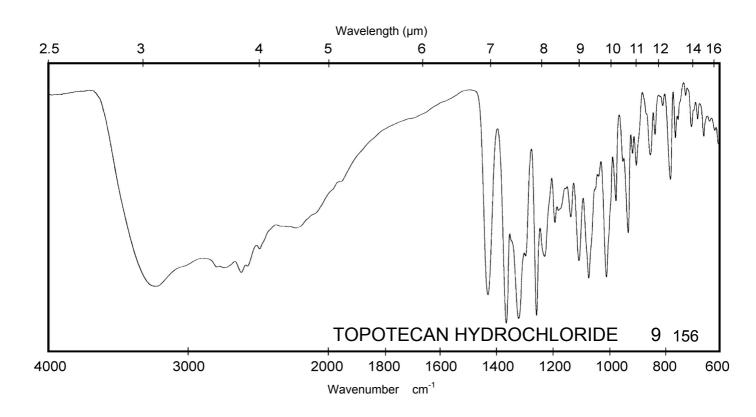




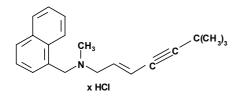


327.9


Concentration 0.4 mg/100 ml 2.6 mg/100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	223 nm 283 nm	222 nm 283 nm	223 nm 283 nm	
E ^{1%} 1cm	2650 242	2410 225	2470 228	
ε	86900 7940	79000 7380	81000 7480	

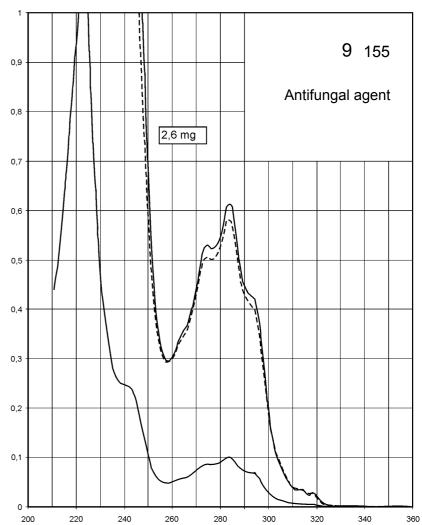

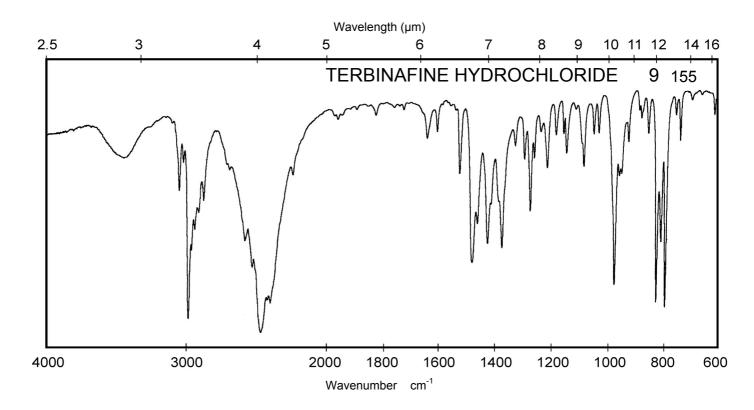




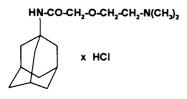


 M_r 316.9


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

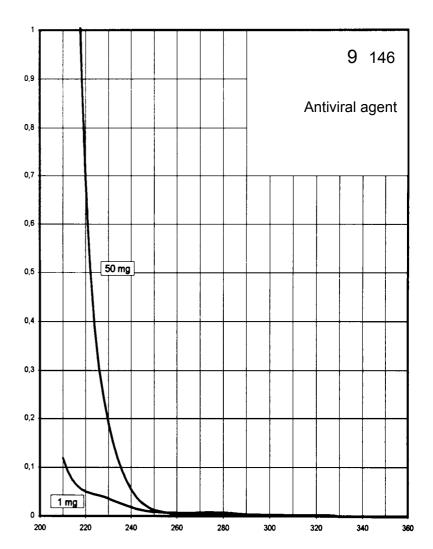

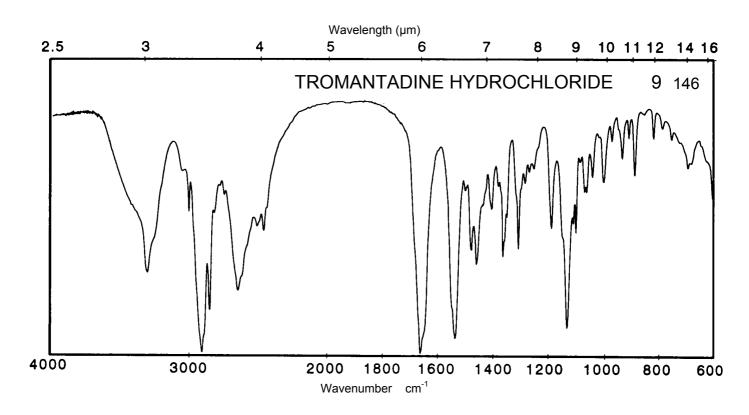



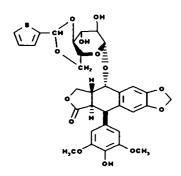



656.7

Concentration


2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm	284 nm	Decom- position observed	288 nm
E ^{1%} _{1cm}	69.3	61.6		124
ε	4500	4000		8200

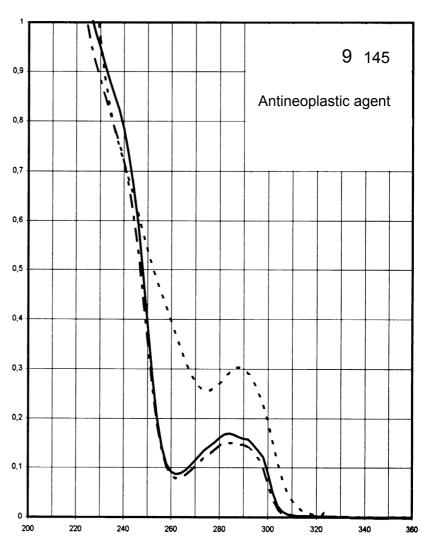

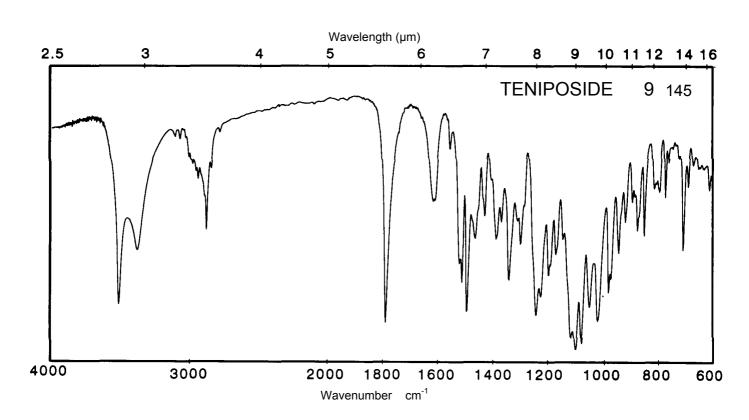





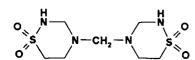


284.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

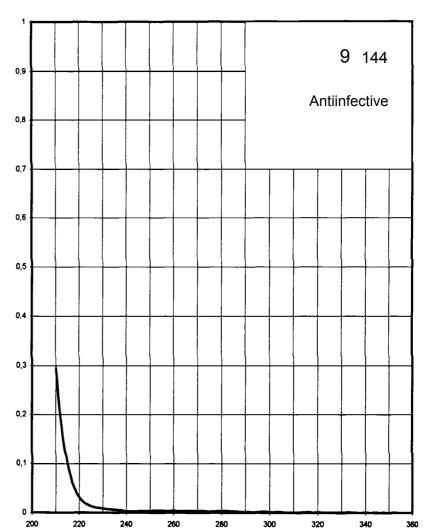

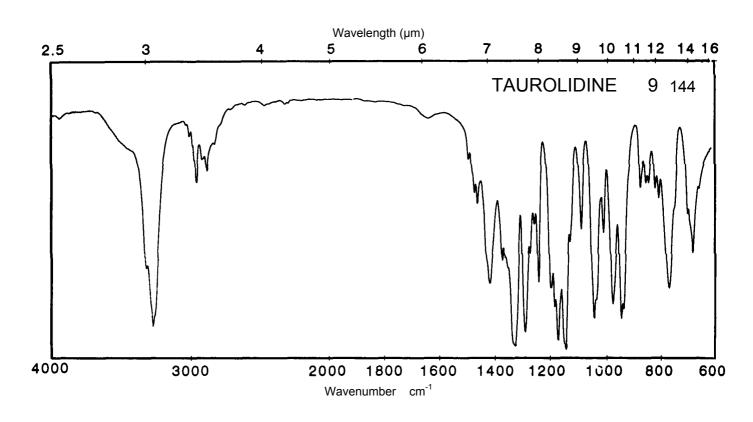





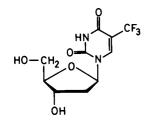


296.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm		261 nm	260 nm
E ^{1%} 1cm	340		345	232
ε	10000		10200	6870

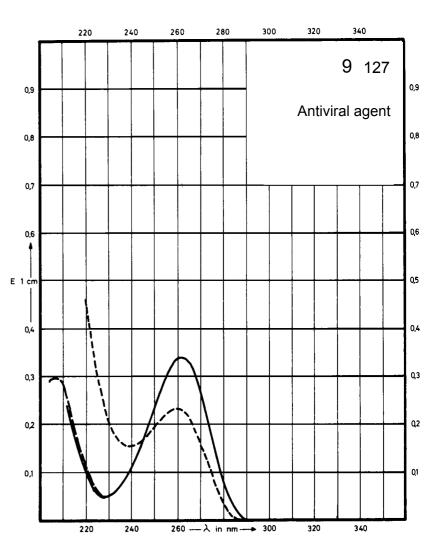

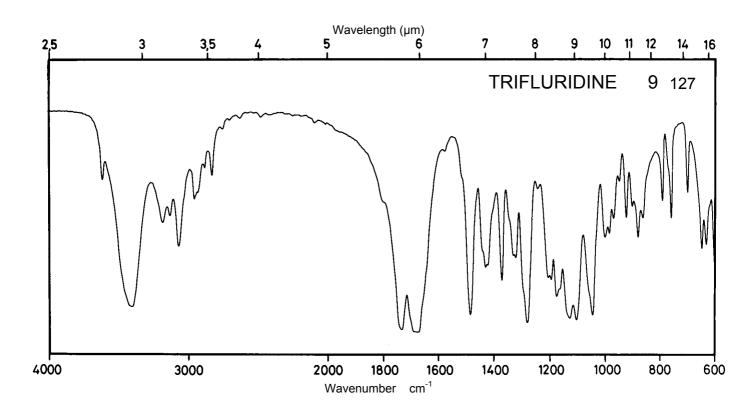





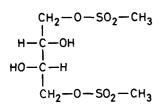


278.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E 1% 1cm				
ε				

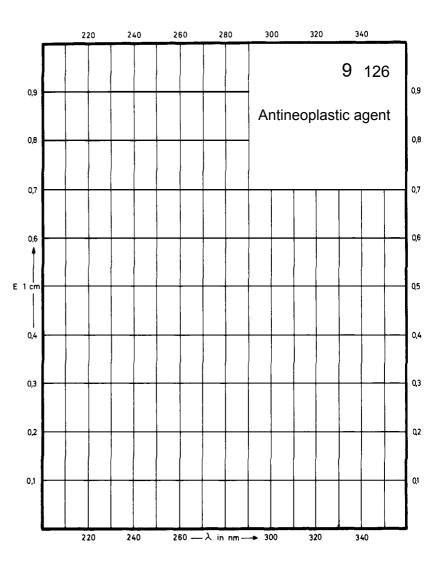

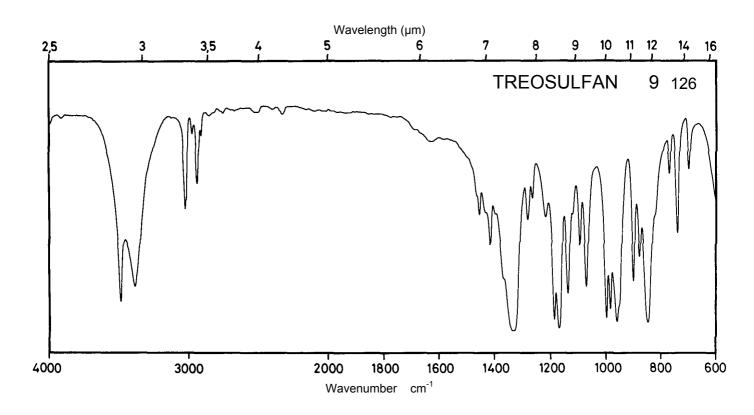





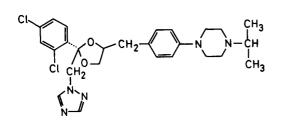

532.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	243 nm 225 nm		225 nm	
E ^{1%} 1cm	268 330		423	
ε	14280 17560		22500	

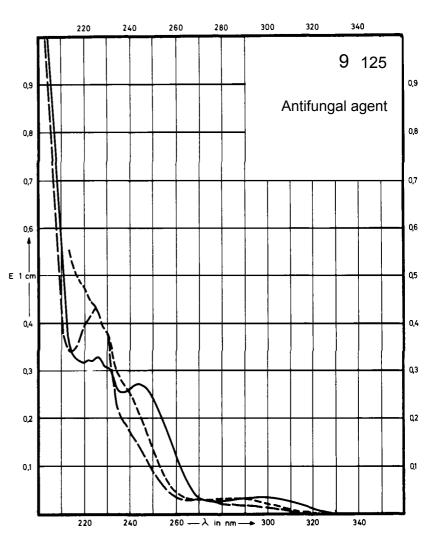

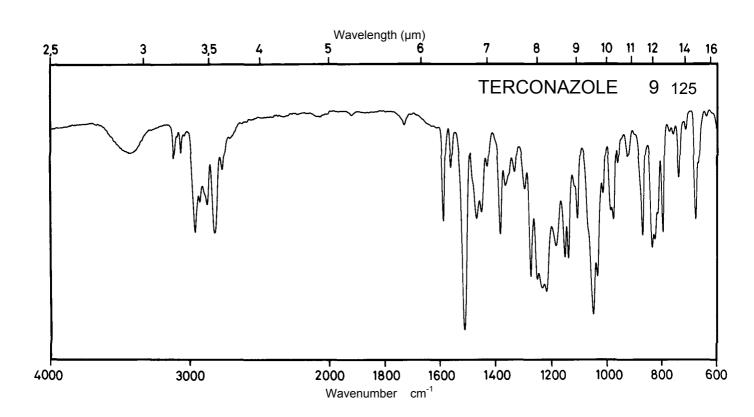


Name TIOCONAZOLE


 M_{r}

387.7


Concentration 1.3 mg / 100 ml 65 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm 271 nm 219 nm		218 nm	
E ^{1%} 1cm	7.2 9.3 454		443	
ε	278 360 17580		17160	

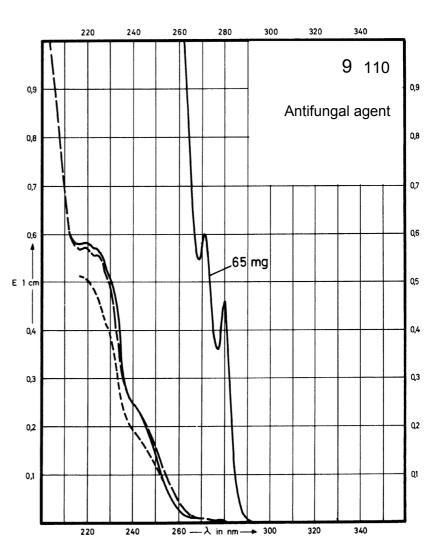

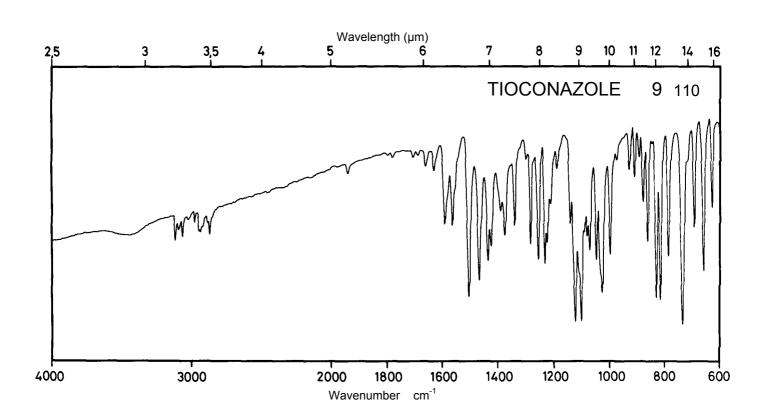



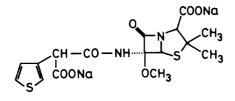



458.4

Concentration


3.1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	(230 nm)		232 nm	
E ^{1%} 1cm	144		152	
3				

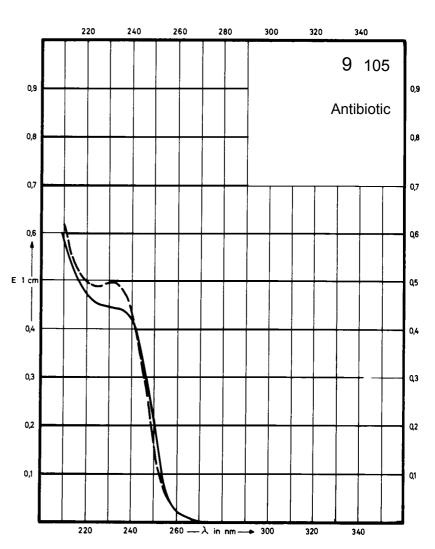

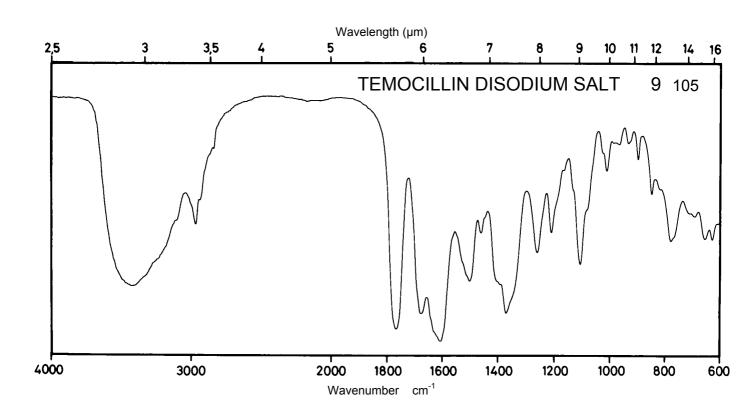





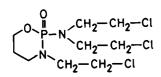


323.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

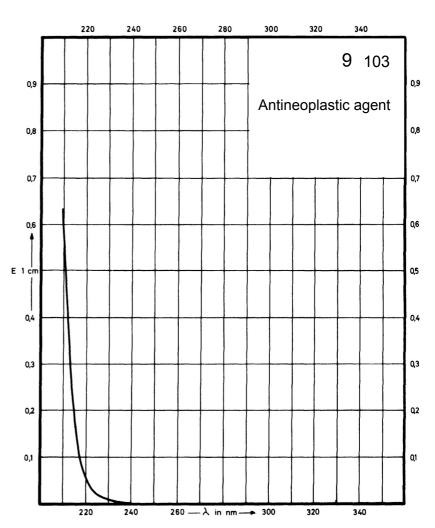

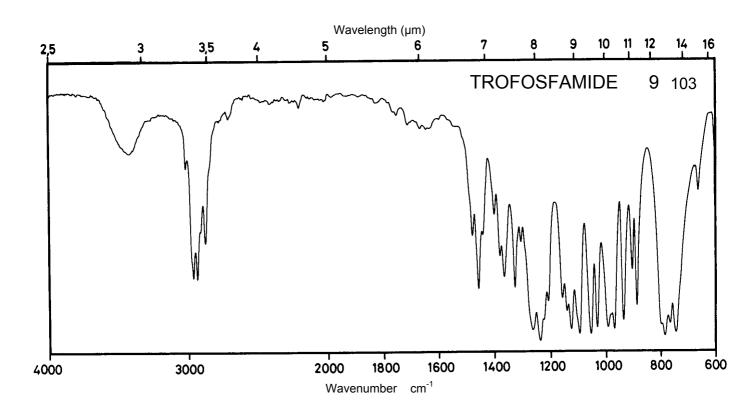





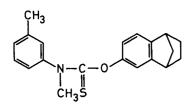


323.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm			
E ^{1%} 1cm	518			
ε	16800			

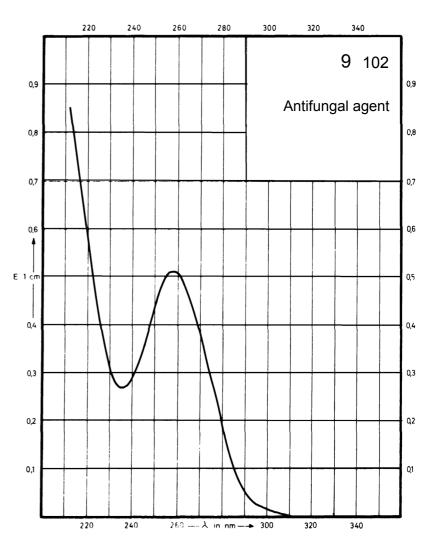

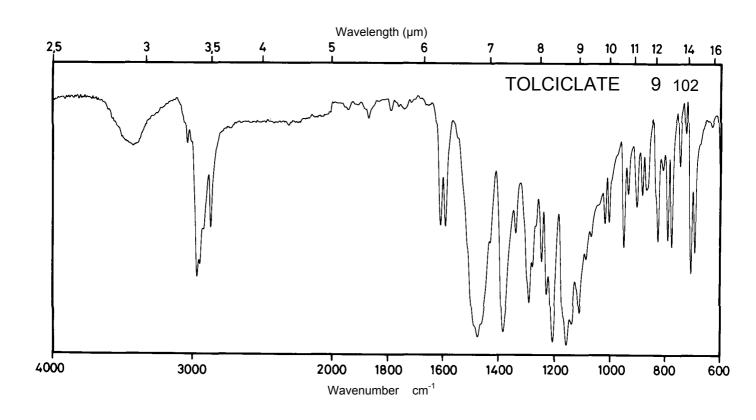





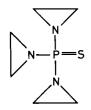


 M_{r} 189.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

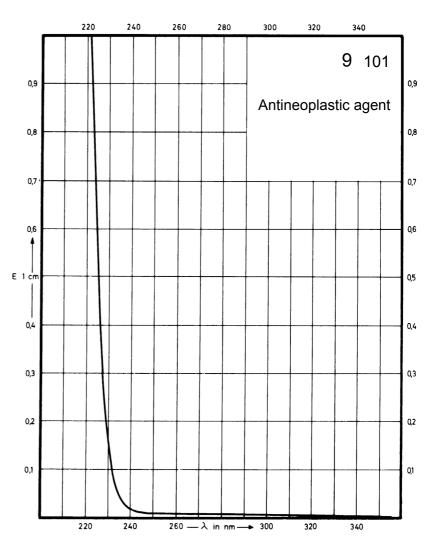

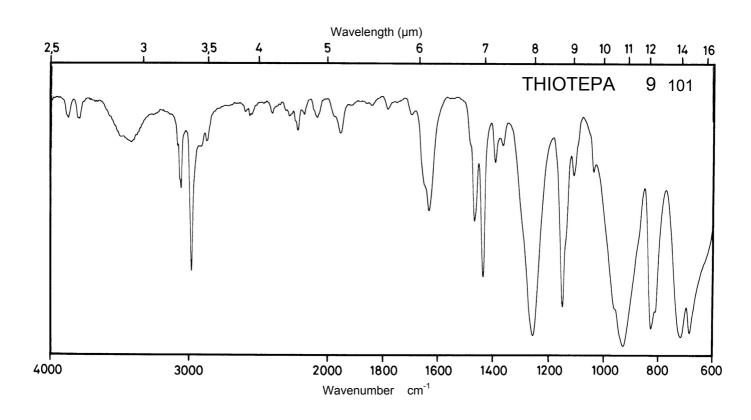




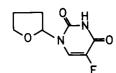


 M_r


200.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm		271 nm	269 nm
E ^{1%} 1cm	434		452	346
ε	8700		9000	6900

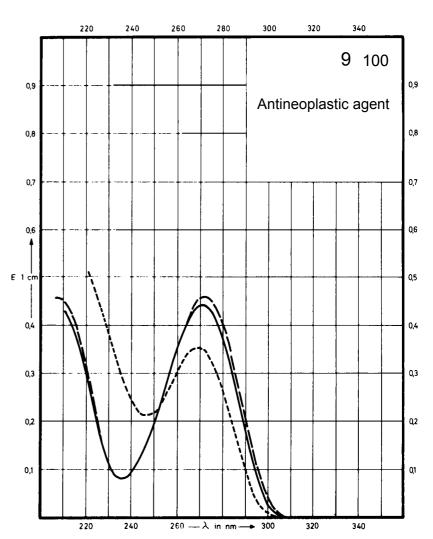

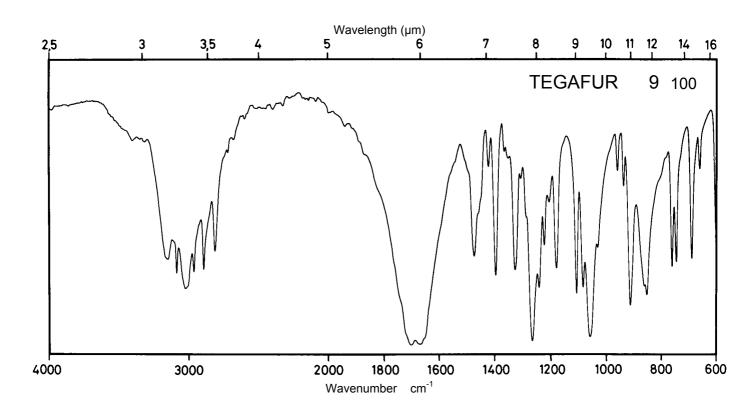





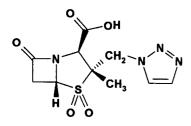


300.3


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				Decom- position observed
E ^{1%} 1cm				
ε				

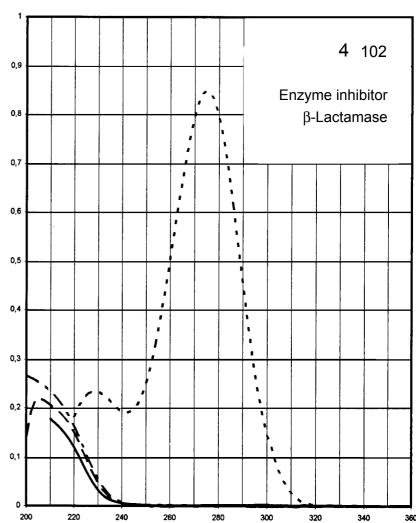


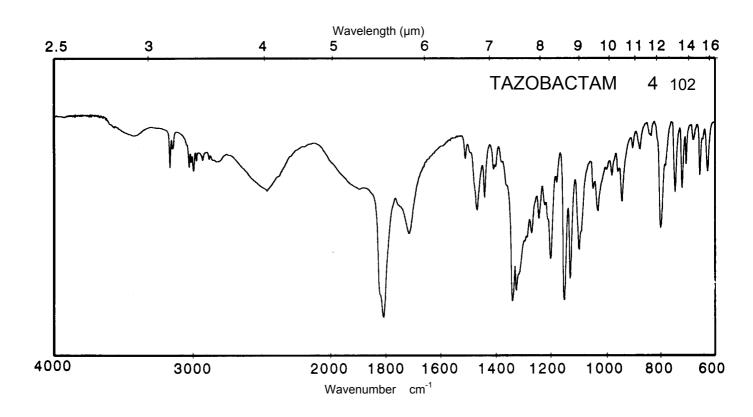


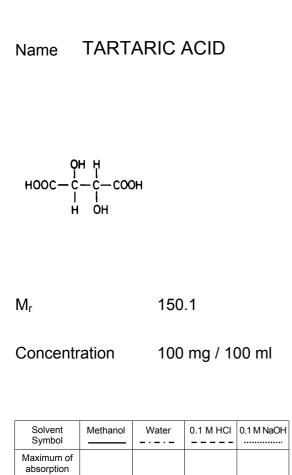


3

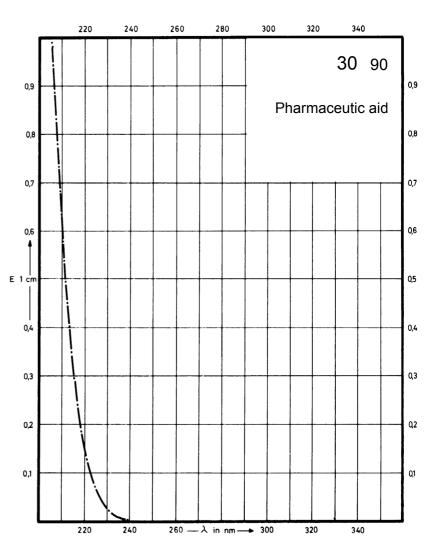


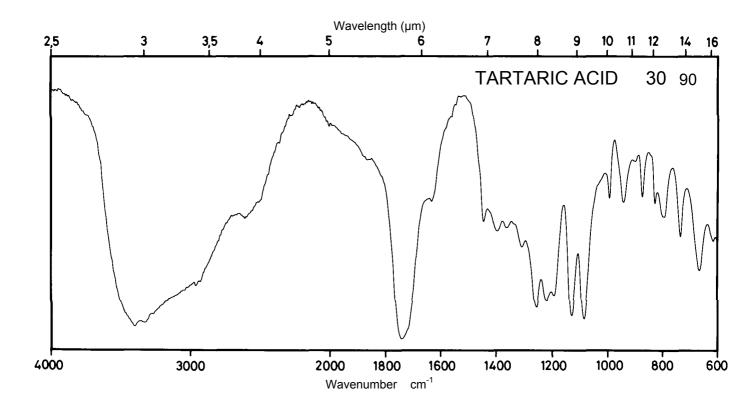

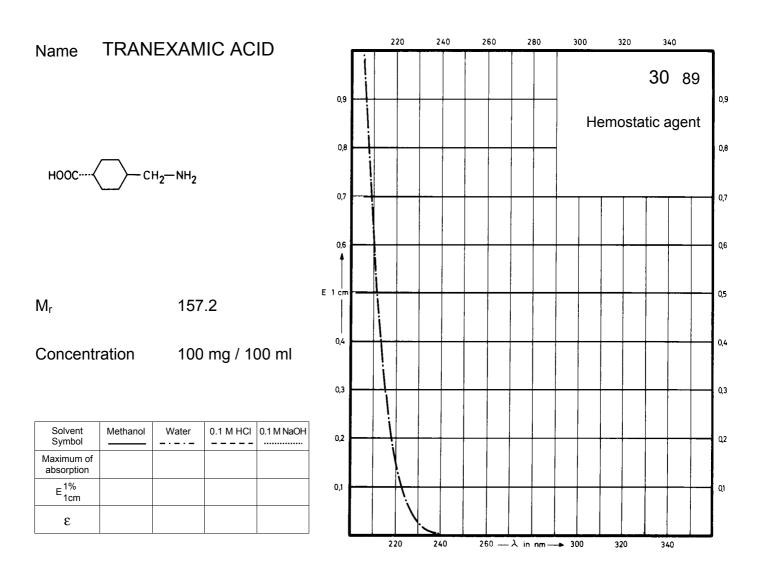

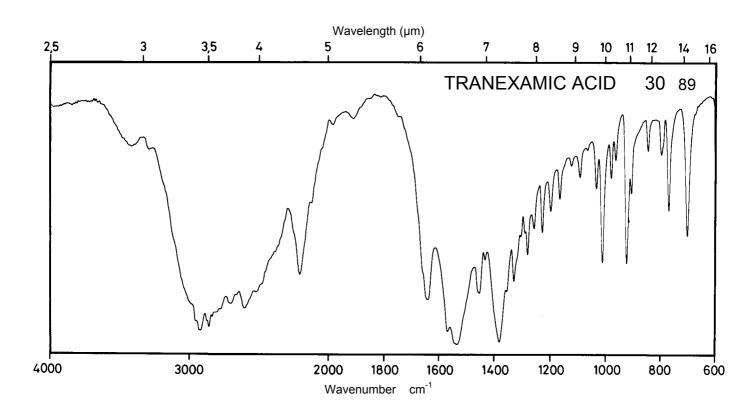


Name TICLOPIDINE HYDROCHLORIDE

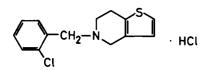


 M_r 300.3


Concentration 2 mg / 100 ml 23 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm 269 nm 233 nm	233 nm	233 nm	
E ^{1%} 1cm	30 34 240	235	235	
ε	890 1000 7200	7100	7100	

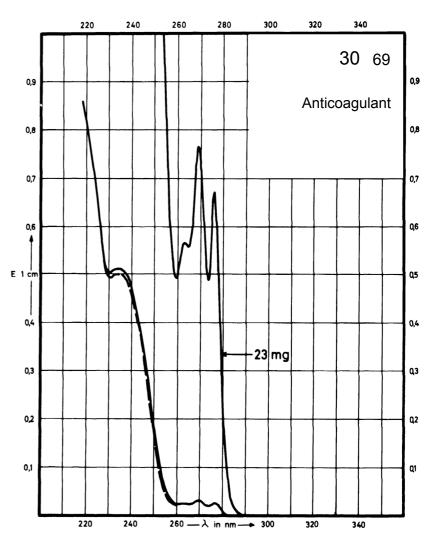

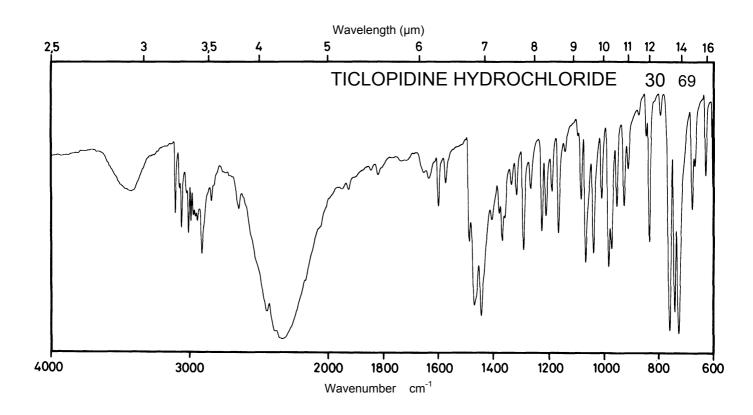





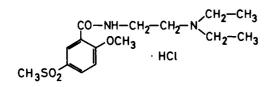


 M_{r} 364.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 235 nm	288 nm	288 nm	285 nm 236 nm
E ^{1%} 1cm	67 383	63	63	62
3	2440 14000	2300	2300	2270

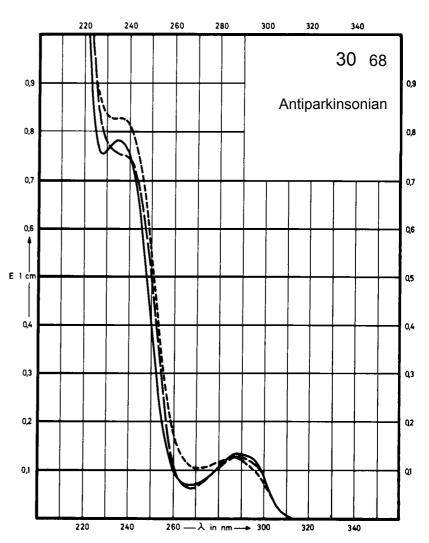

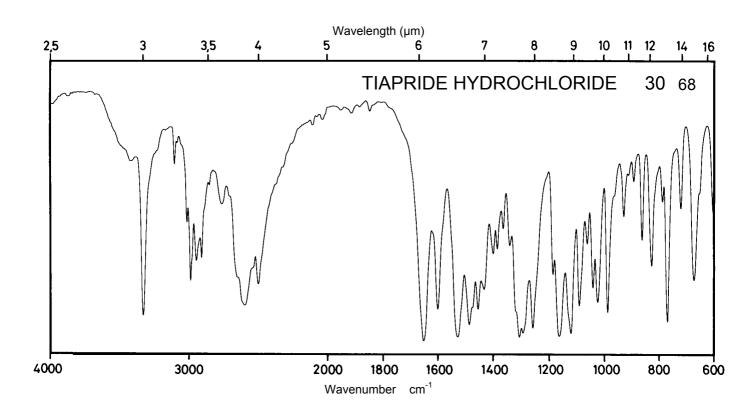





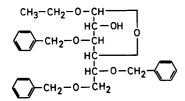


478.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm 253 nm			
E ^{1%} 1cm	13.0 12.0			
ε	620 570			

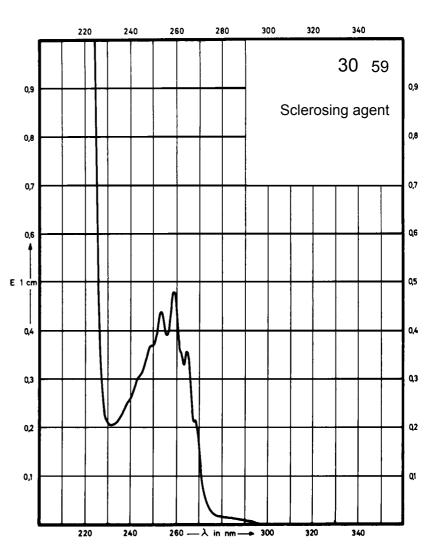

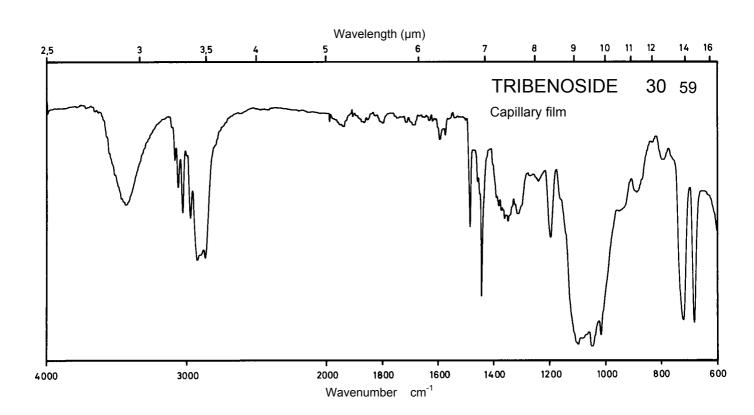





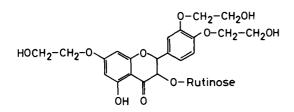


742.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	352 nm 255 nm		350 nm 254 nm	370 nm 280 nm
E ^{1%} 1cm	257 315		244 300	114 314
ε	19090 23400		18110 22300	8440 23320

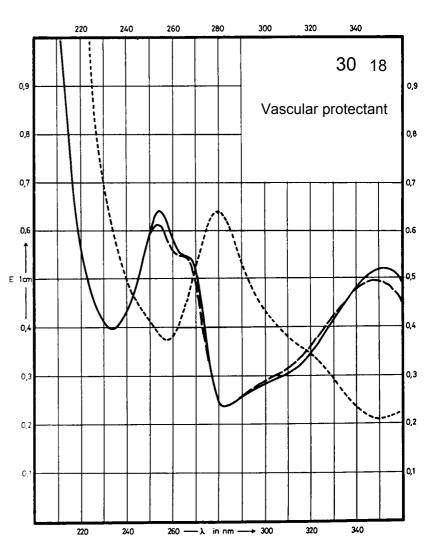

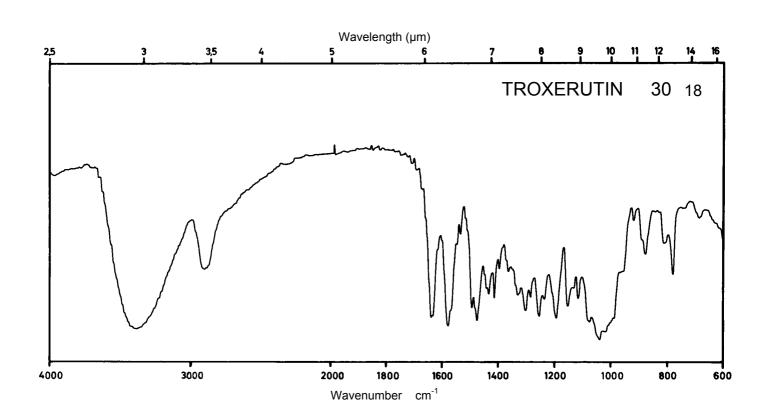




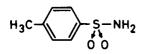


 M_r 171.2


Concentration 1 mg / 100 ml 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm 224 nm		261 nm 224 nm	260 nm
E ^{1%} 1cm	23 725		31 757	27
3	390 12410		530 12960	460

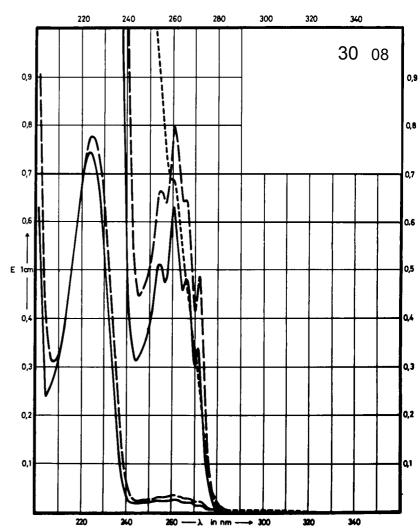

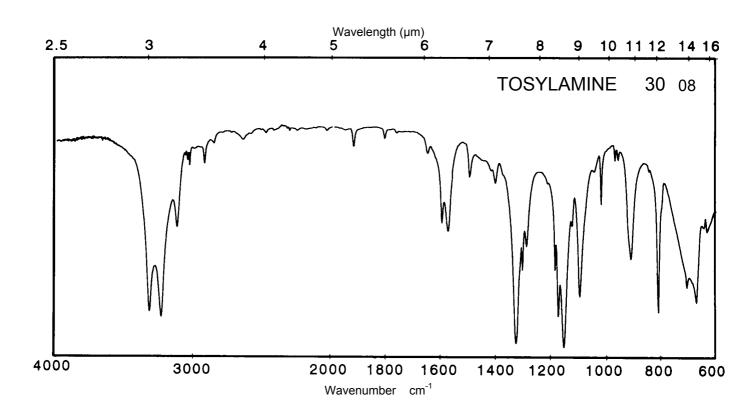





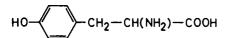


181.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			274 nm	294 nm
E ^{1%} 1cm			76	130
ε			1380	2360

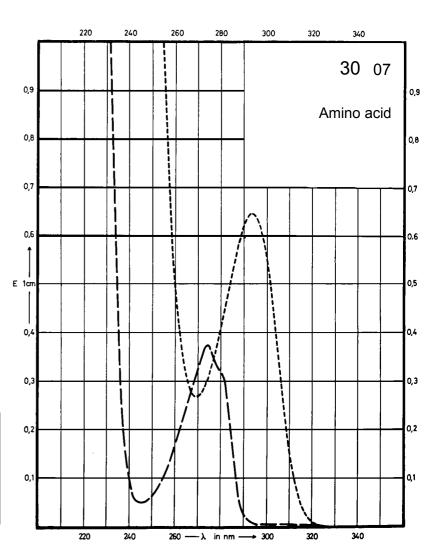

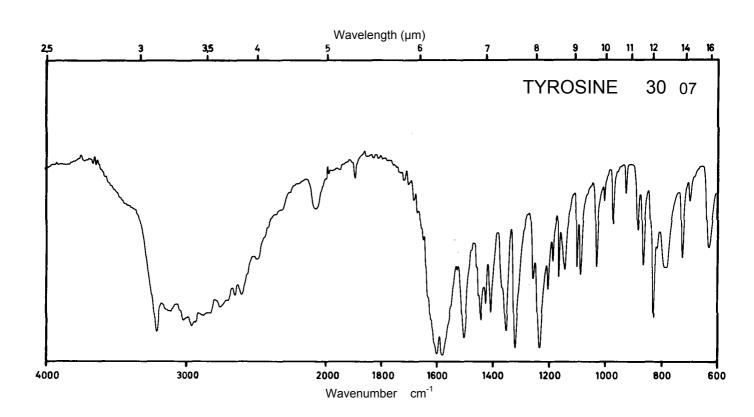




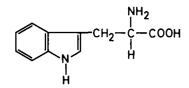



204.2

Concentration

1 mg / 100 ml

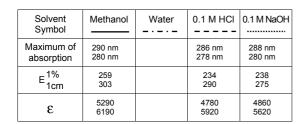

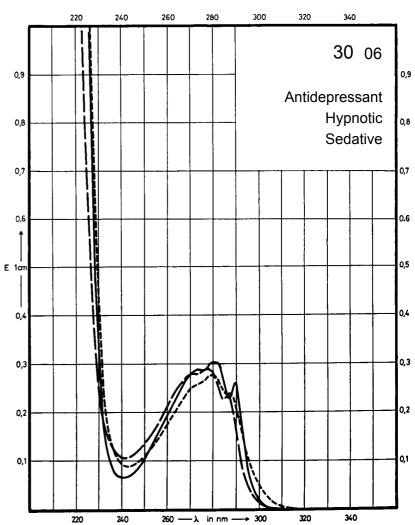

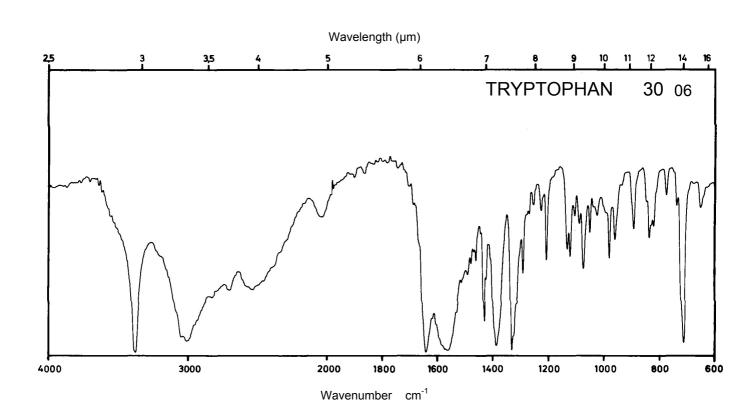


Name TOLDIMFOS SODIUM SALT


 M_{r}

221.2


Concentration 0.5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		275 nm 267 nm	268 nm
E 1% 1cm	1085		48 48	770
ε	24000		1060 1060	17030

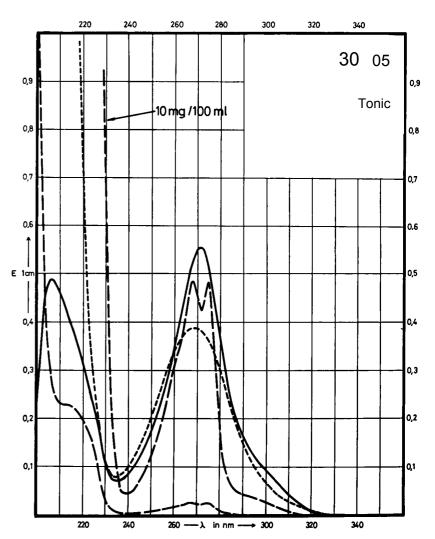

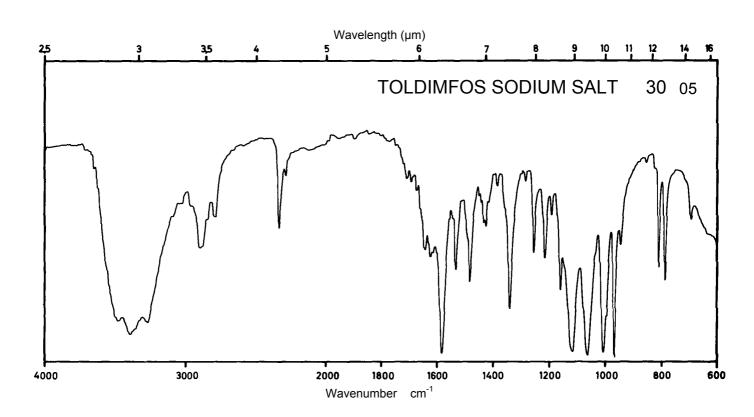


Name TRYPTAMINE HYDROCHLORIDE

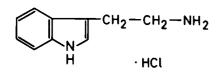


 M_{r}

196.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm		278 nm	288 nm
E ^{1%} 1cm	311		288	281
ε	6120		5660	5530

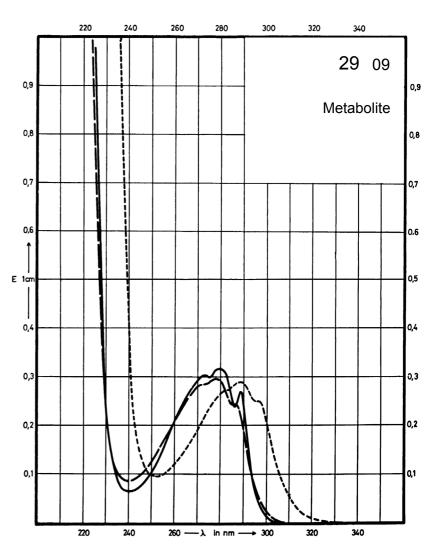

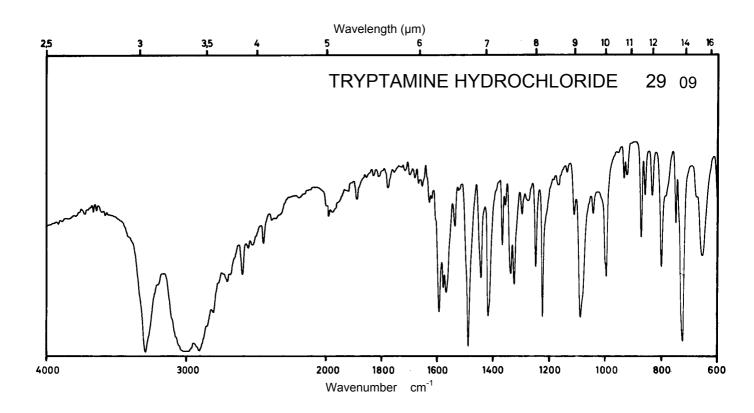



Name TOCOPHEROL ACETATE, α -

 M_{r}


472.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm			
E ^{1%} 1cm	43			
3	2030			

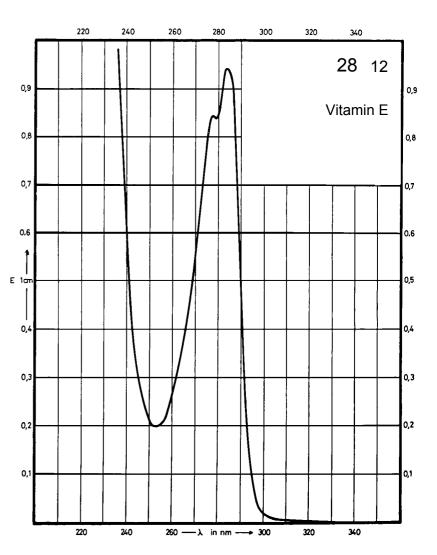

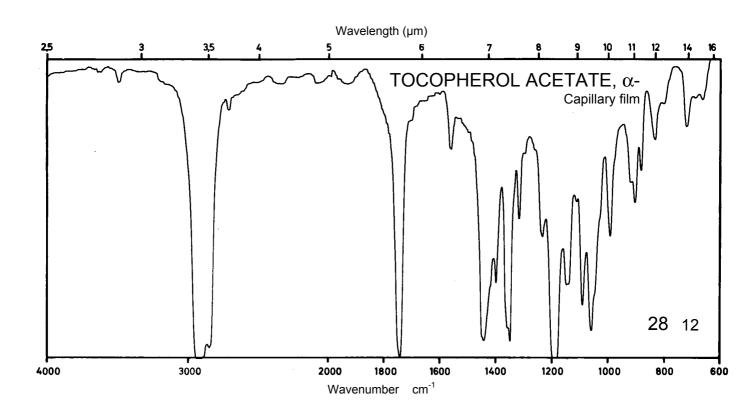




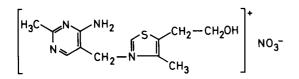


 M_r

Concentration 1.5 mg / 100 ml

327.4

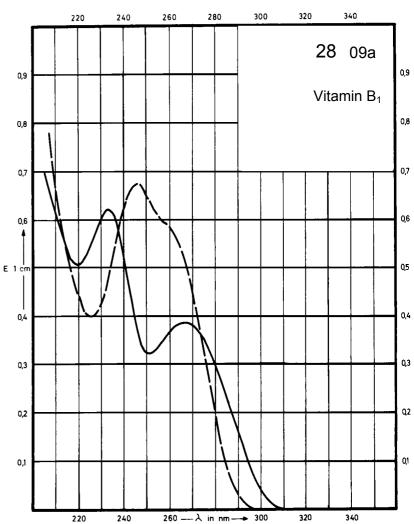

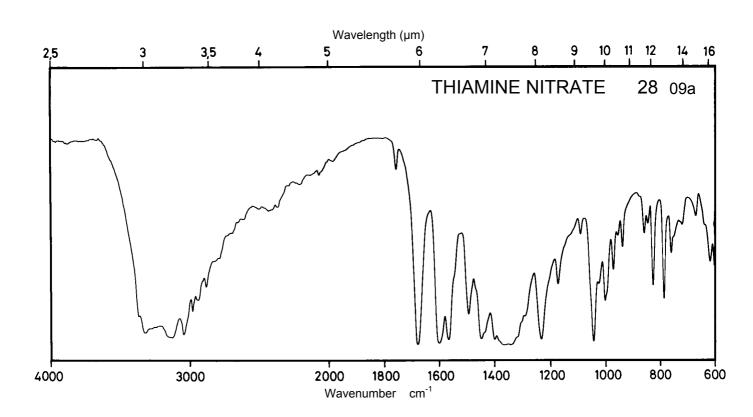



Name THIAMINE HYDROCHLORIDE

 M_{r}


337.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm 238 nm		246 nm	232 nm
E ^{1%} 1cm	286 352		413	505
3	9650 11870		13930	17030

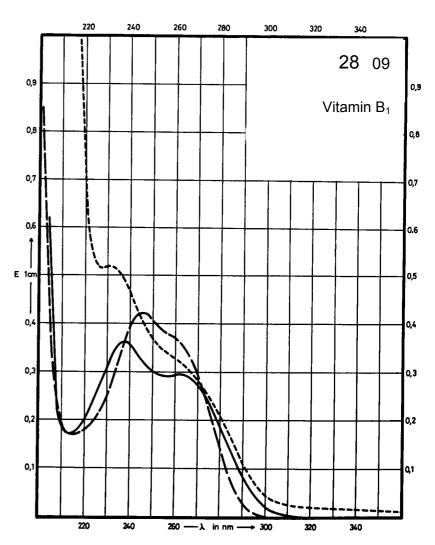

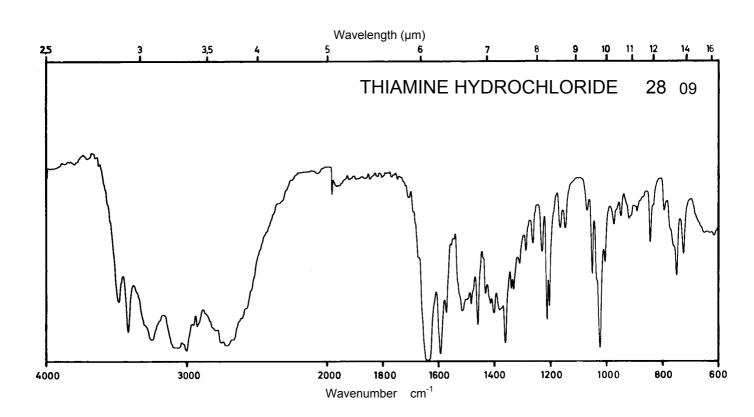

Name TRAMAZOLINE HYDROCHLORIDE

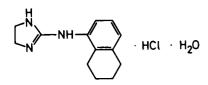



 M_{r}

269.8


Concentration 5 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm 265 nm		272 nm 265 nm	
E ^{1%} _{1cm}	21.7 25.3		25.0 28.5	
ε	590 680		680 770	

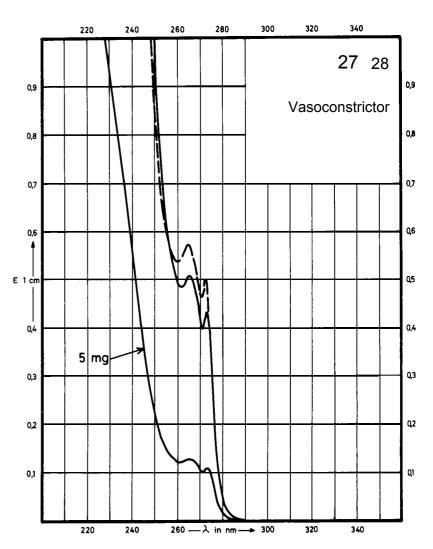

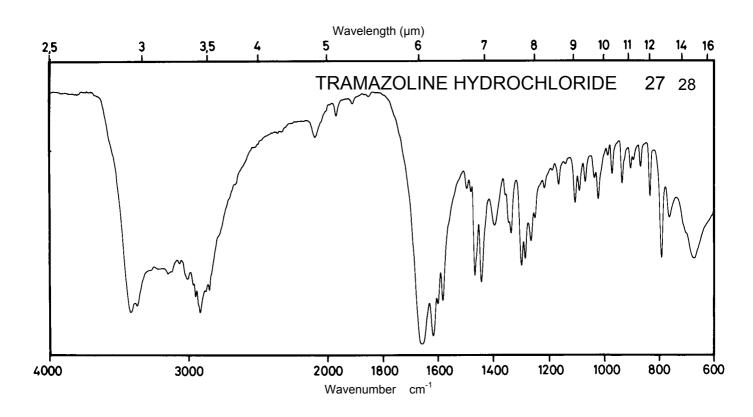


Name TETRYZOLINE HYDROCHLORIDE

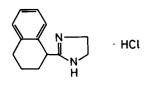


 M_r 200.3

Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm	272 nm	272 nm	273 nm
	266 nm	265 nm	265 nm	266 nm
E ^{1%} 1cm	15.4	15.8	16.0	16.0
	16.0	17.3	17.4	18.0
ε	309	317	320	320
	320	346	348	360

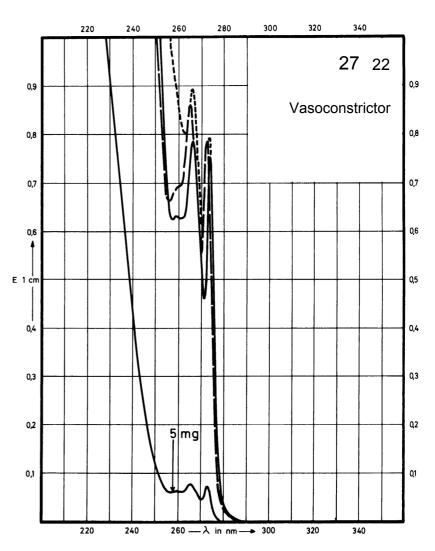

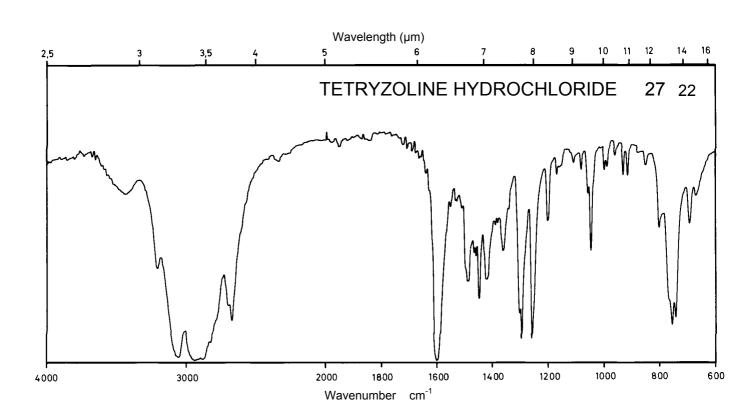



Name TYRAMINE HYDROCHLORIDE

 M_{r}


173.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm		273 nm	294 nm
E ^{1%} 1cm	100		86	144
ε	1740		1490	2500

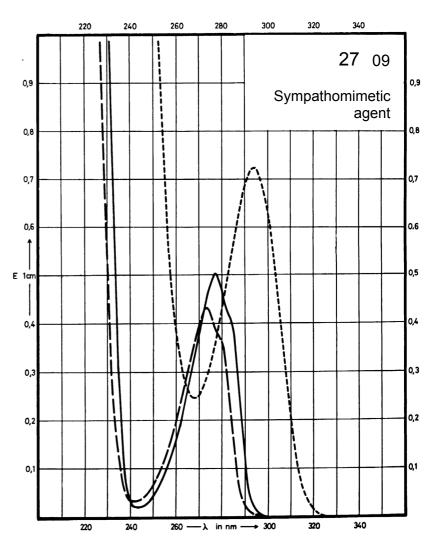

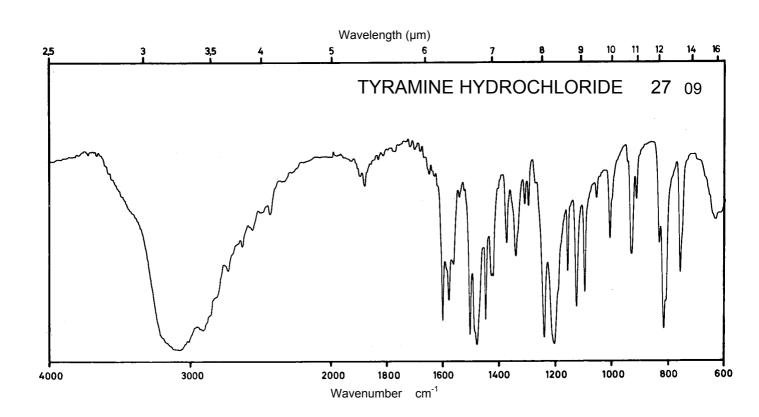




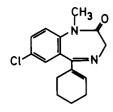


 M_r 288.8


Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	307 nm 226 nm	227 nm	345 nm 284 nm 239 nm	305 nm 227 nm
E ^{1%} 1cm	76 1020	1060	149 452 940	84 1110
ε	2190 29460	30600	4300 13050 27150	2430 32060

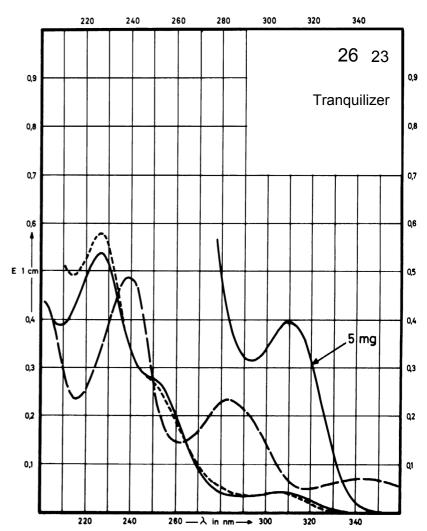

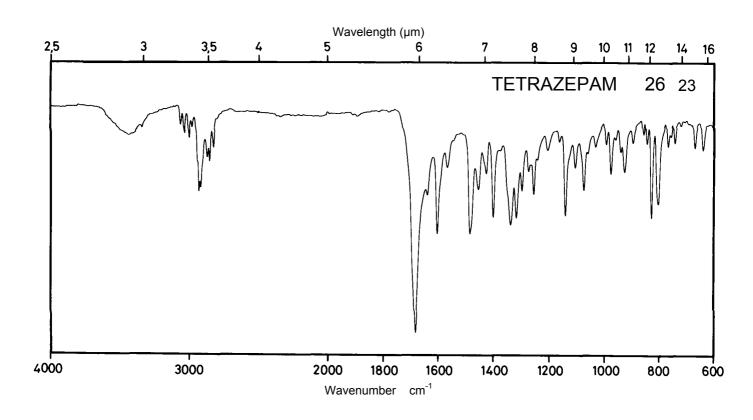




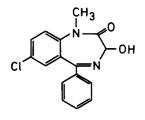


 M_r 300.7

Concentration 0.5 mg / 100 ml 1 mg / 100 ml 5.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	314 nm 230 nm		358 nm 284 nm 234 nm	310 nm
E ^{1%} 1cm	76 1080		68 283 980	81
3	2300 32480		2050 8500 29420	2440

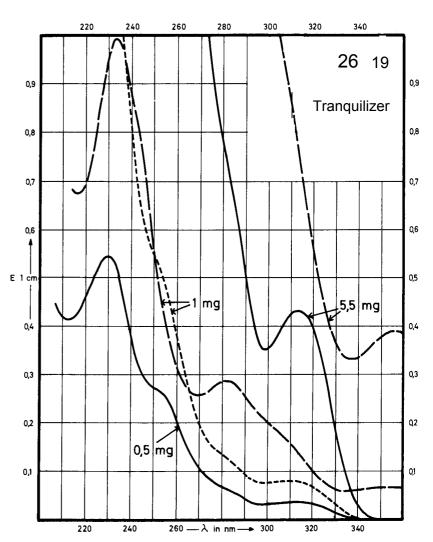

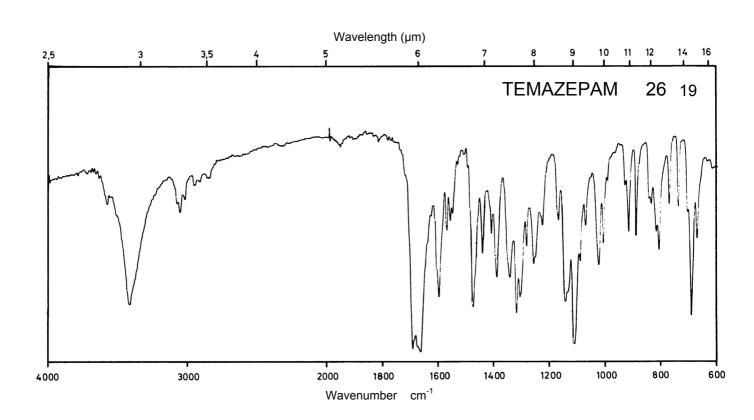

Name TRIMETOZINE


$$H_3C-0$$
 H_3C-0
 $CO-N$


 M_{r}

281.3


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm			
E ^{1%} 1cm	259			
3	7290			

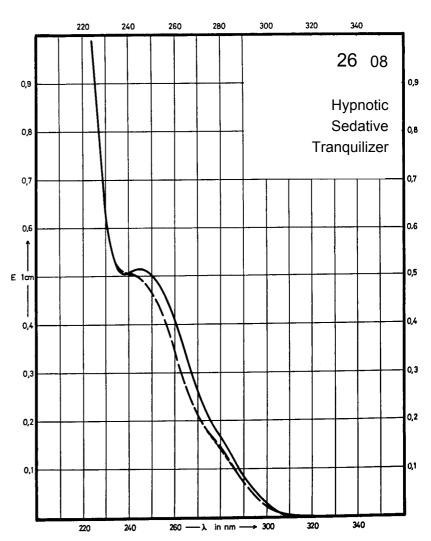

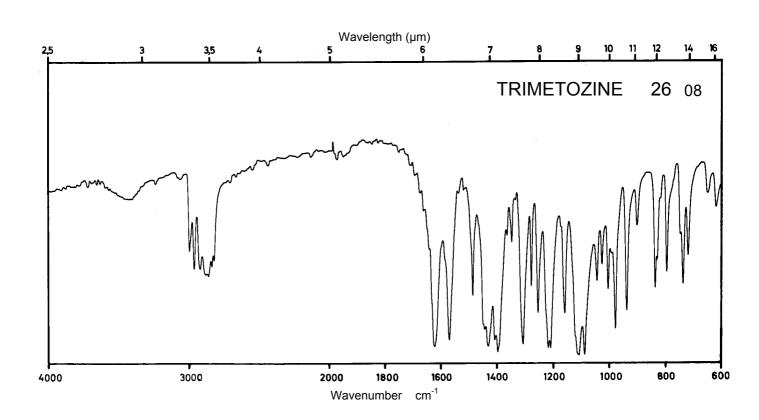




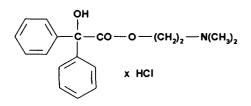



335.8


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm 252 nm			
E ^{1%} 1cm	13.8 16.5 14.8			
ε	460 550			

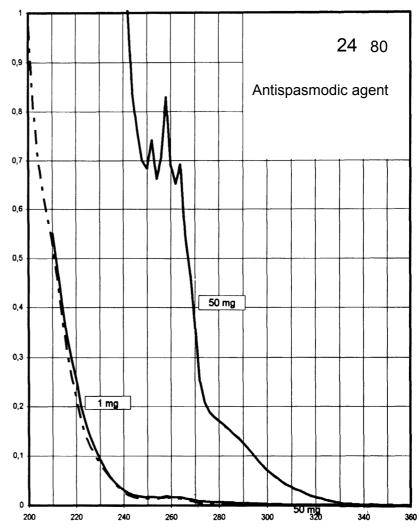


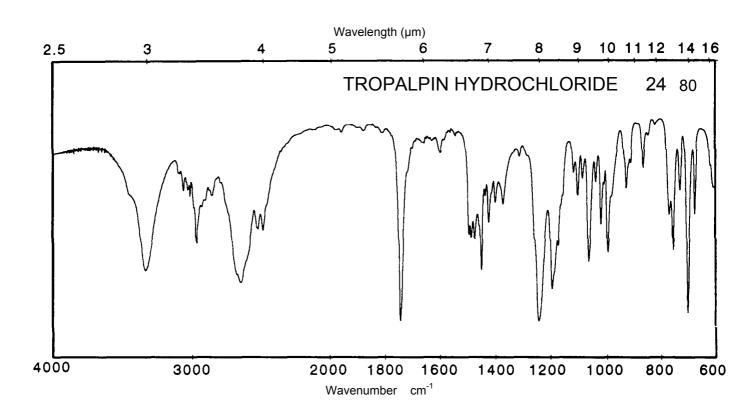

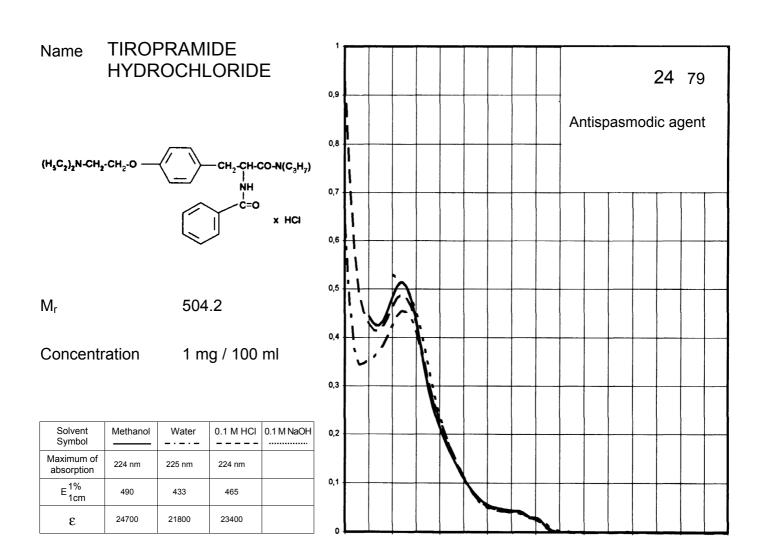

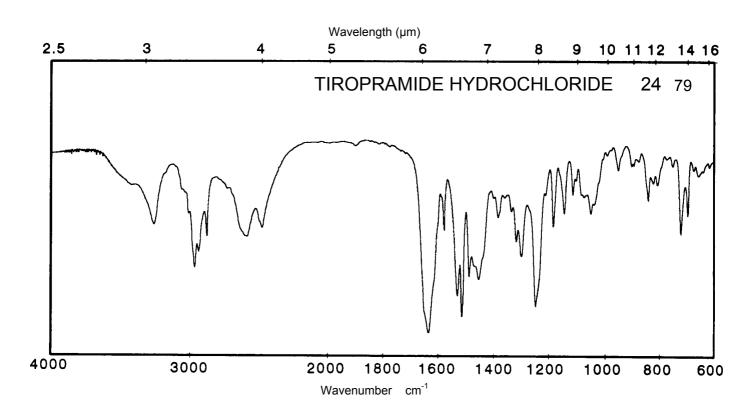

Name TULOBUTEROL HYDROCHLORIDE

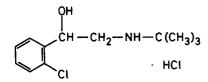



 M_{r}

264.2


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	261 nm		265 nm 261 nm	
E ^{1%} 1cm	13.1		14.9 14.8	
ε	345		393 392	

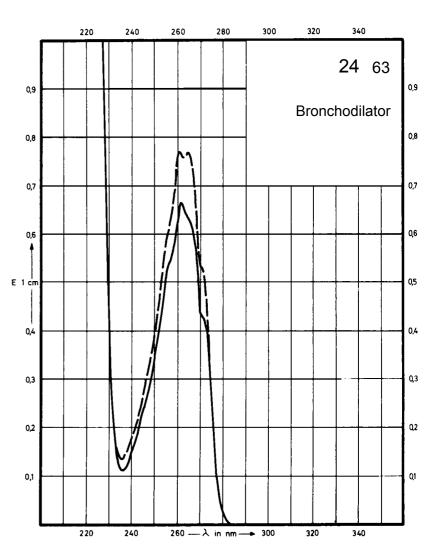

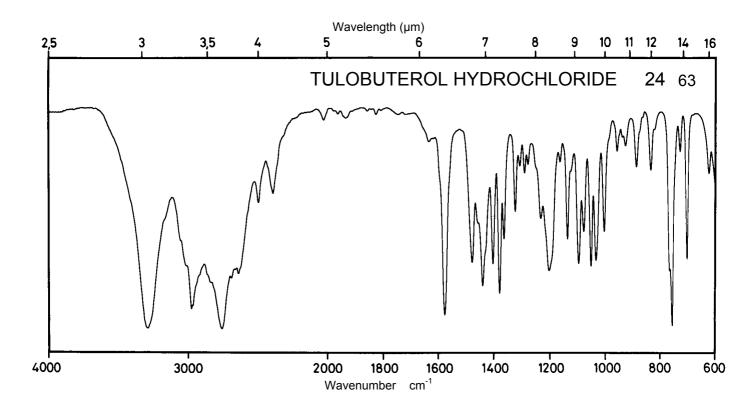

Name TIZANIDINE HYDROCHLORIDE

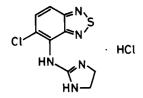



 M_{r}

290.2


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			320 nm 227 nm	319 nm 238 nm
E ^{1%} 1cm			536 666	355 822
ε			15550 19300	10300 23900

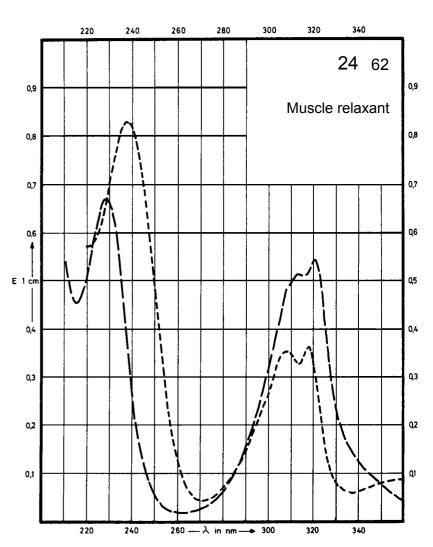

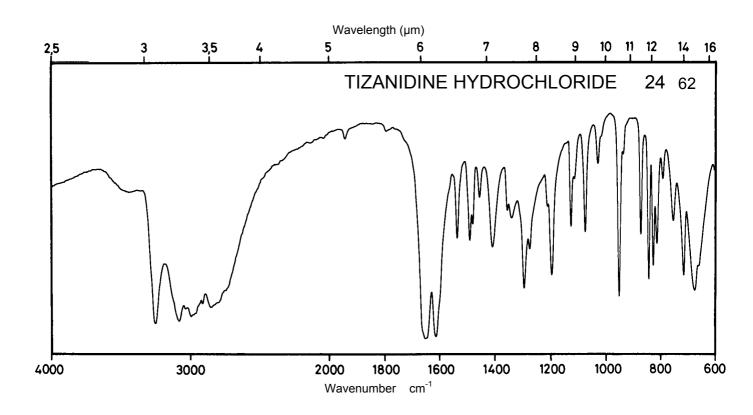

Name THIOCOLCHICOSIDE

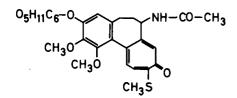



 M_{r}

563.6


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	368 nm 288 nm 256 nm	377 nm 259 nm	375 nm 259 nm	378 nm 259 nm
E ^{1%} 1cm	327 220 376	330 394	330 393	333 398
ε	18410 12430 21200	18600 22200	18600 22160	18770 22430

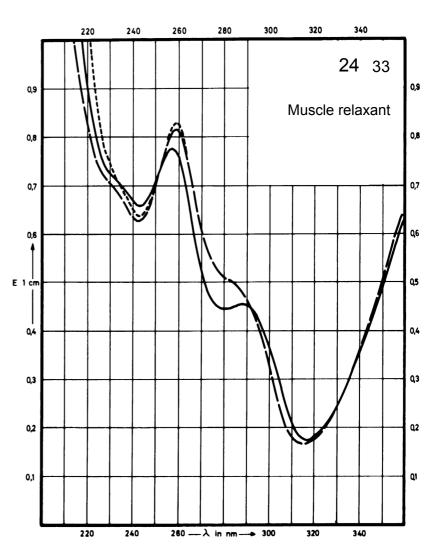

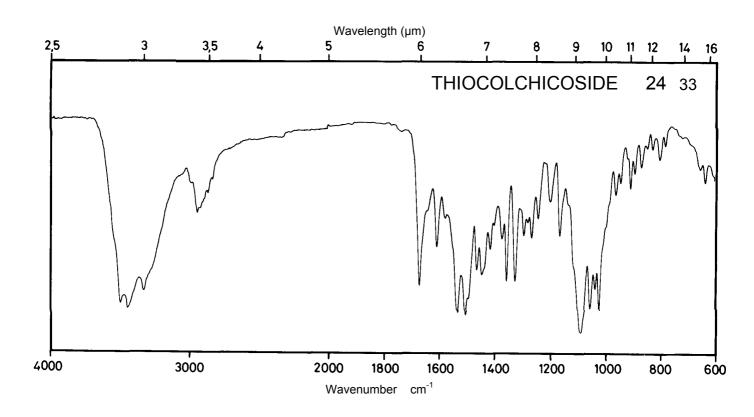




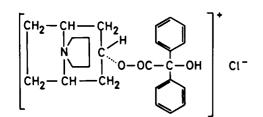


428.0


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm 259 nm 253 nm		265 nm 259 nm 253 nm	
E ^{1%} 1cm	8.07 10.3 8.83		8.53 10.8 8.97	
ε	346 442 378		365 461 384	

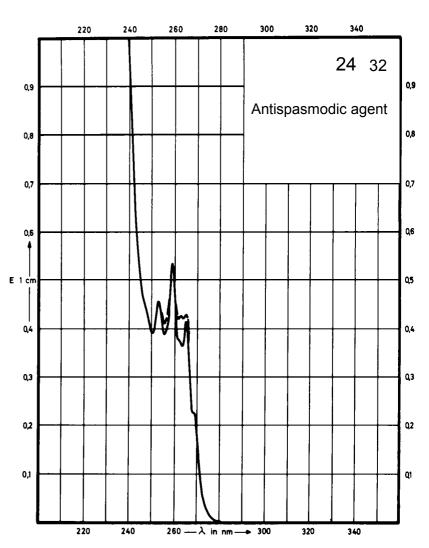

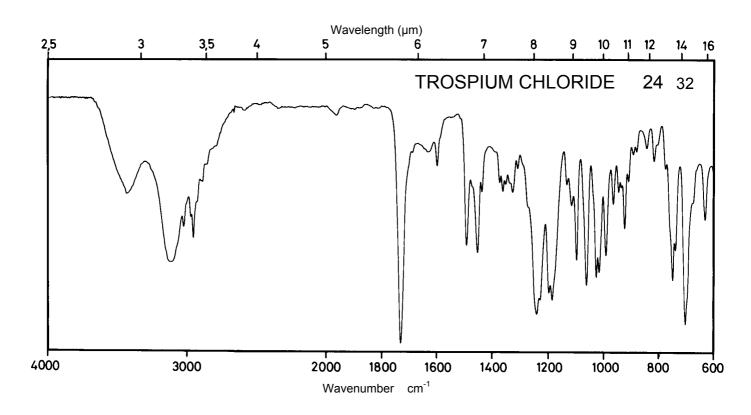




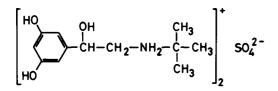



 M_r 548.7


Concentration 2 mg / 100 ml 4 mg / 100 ml 8 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	283 nm 278 nm	276 nm	276 nm	296 nm
E ^{1%} 1cm	75 75	68	68	123
3	4120 4120	3730	3730	6750

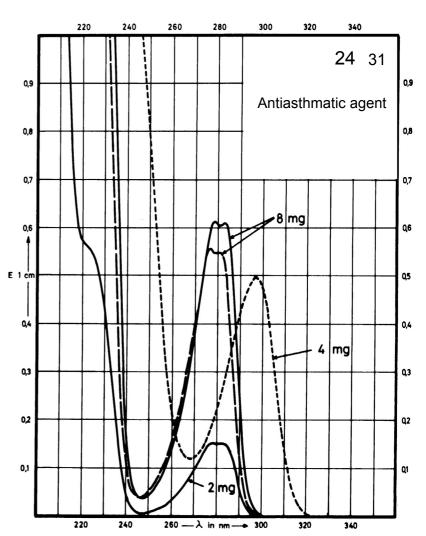


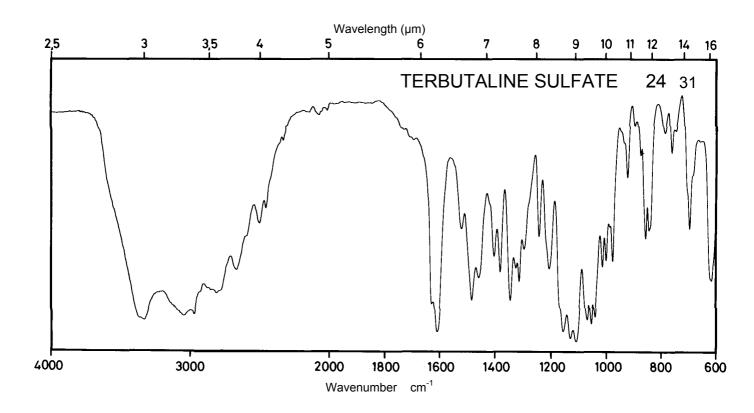

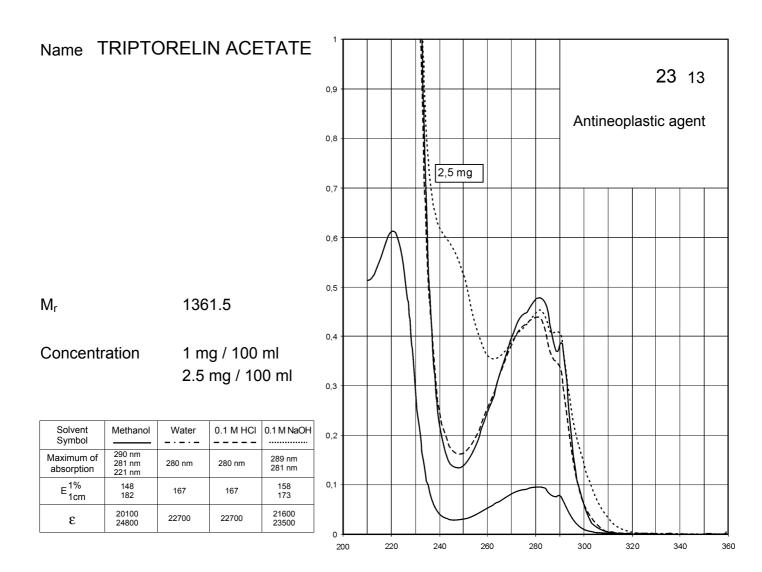


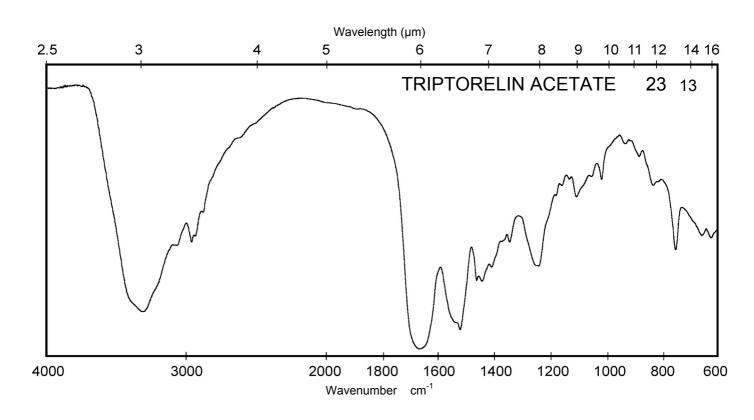

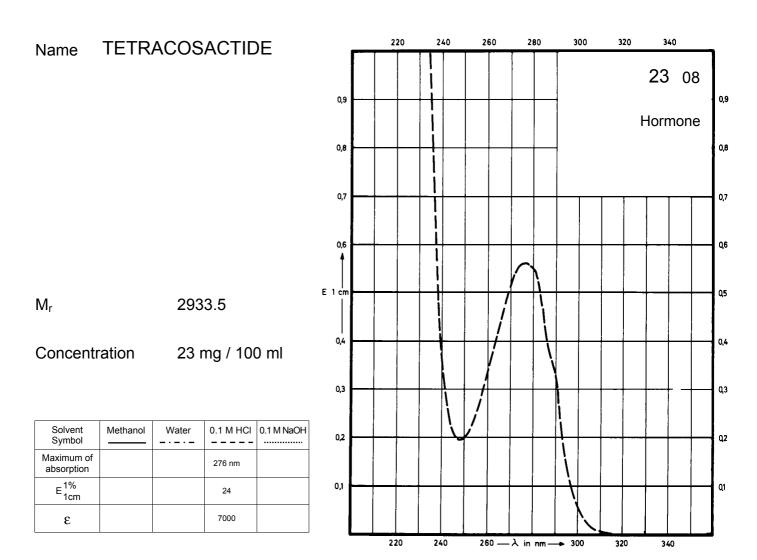

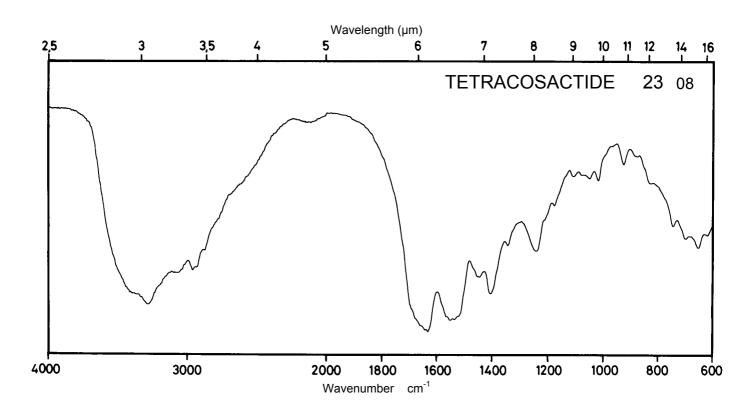




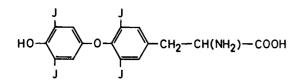


776.9


Concentration 9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				324 nm
E ^{1%} 1cm				75
3				5830

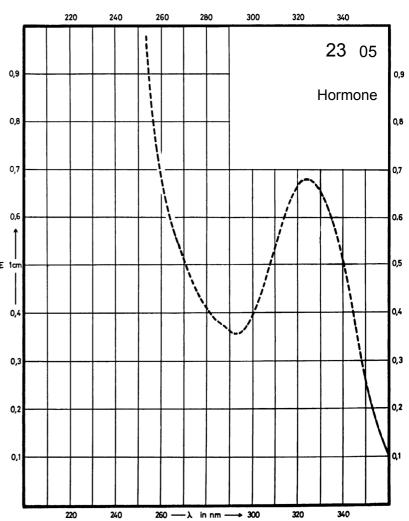

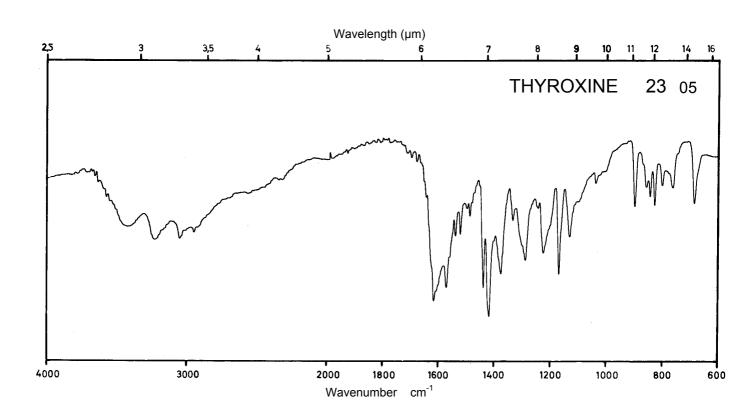



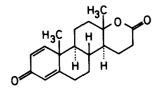



300.4

Concentration


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm	246 nm		
E ^{1%} 1cm	545	558		
3	16370	16760		

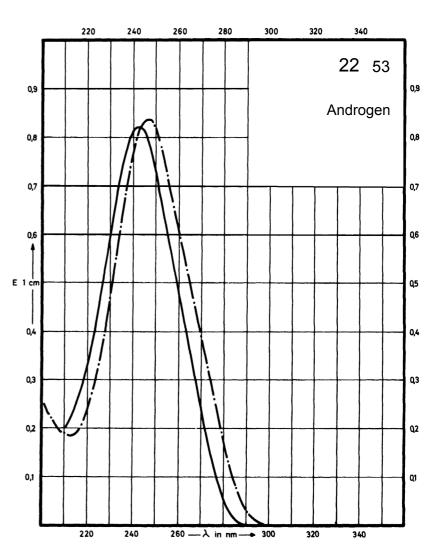

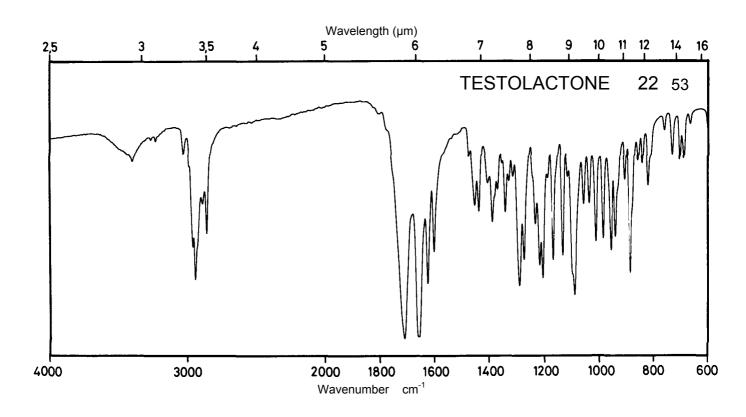


Name TESTOSTERONE

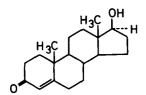


 M_r 288.4


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	241 nm		249 nm	249 nm
E ^{1%} 1cm	564		560	560
3	16270		16150	16150

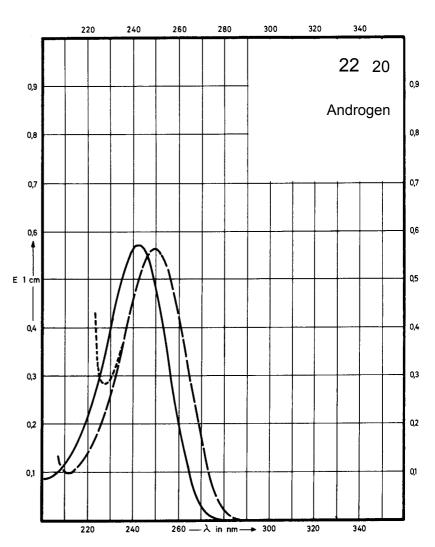

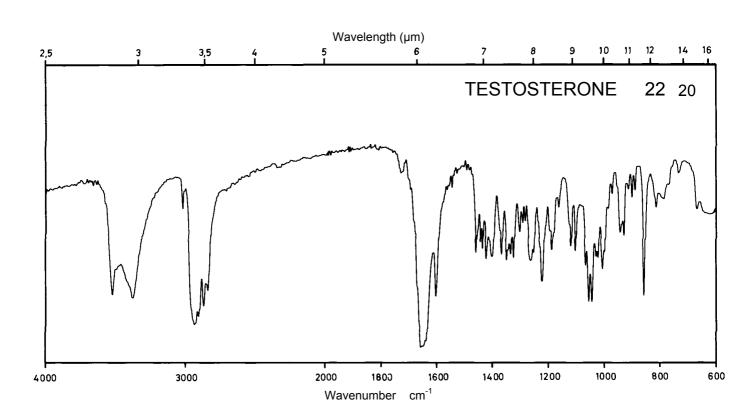


Name TESTOSTERONE 17β -CIPIONATE

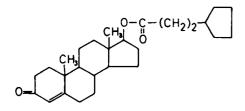


 M_r 413.6

Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm			
E ^{1%} 1cm	407			
3	16830			

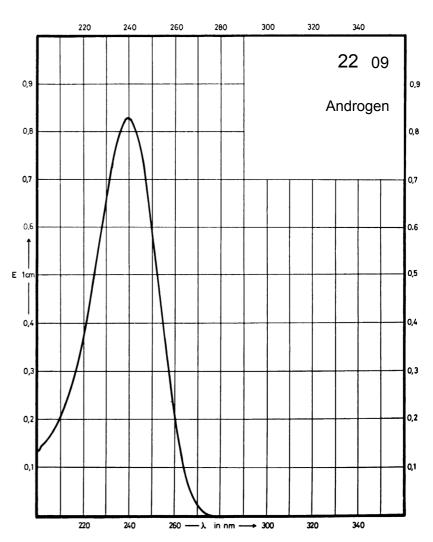

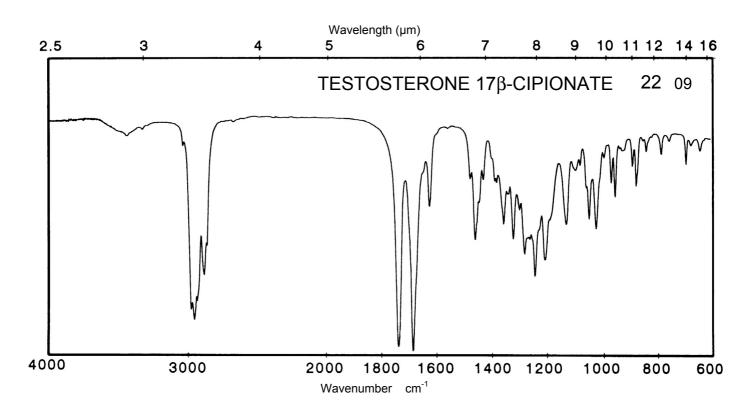


Name THIOPROPERAZINE MESILATE


$$\begin{array}{c|c} S & CH_3 \\ N & SO_2 - N - CH_3 \\ CH_2 - CH_2 - CH_2 - N & N - CH_3 \\ & 2 \ CH_3SO_3H \end{array}$$

 M_r 638.8

Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	318 nm 265 nm		314 nm 263 nm	314 nm 265 nm
E ^{1%} 1cm	55 570		48 480	52 500
ε	3500 36400		3100 30700	3300 31900

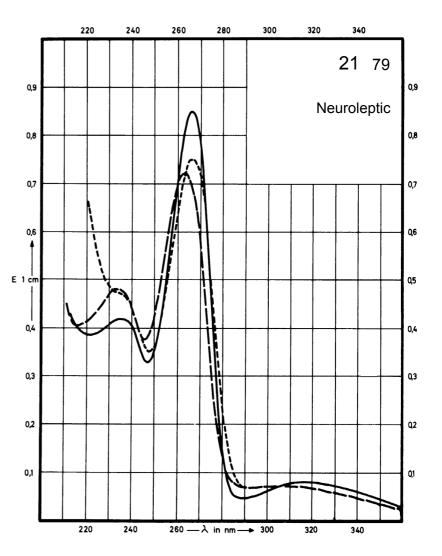

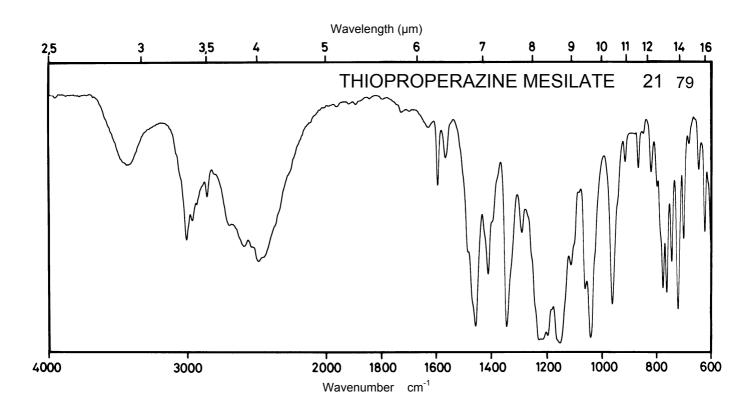


Name TRIFLUPROMAZINE HYDROCHLORIDE


 M_r

292.8

Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	308 nm 258 nm		305 nm 256 nm	
E ^{1%} 1cm	91 853		91 783	
3	2680 24980		2680 22920	

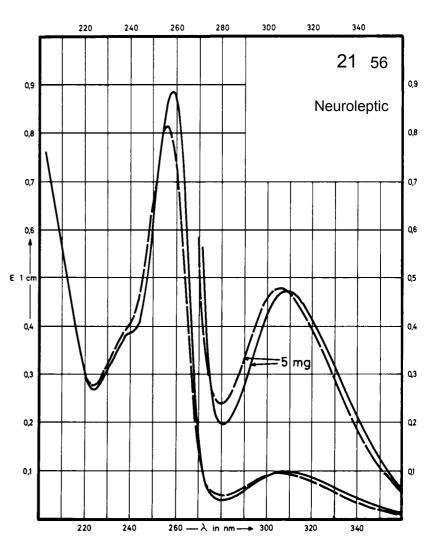

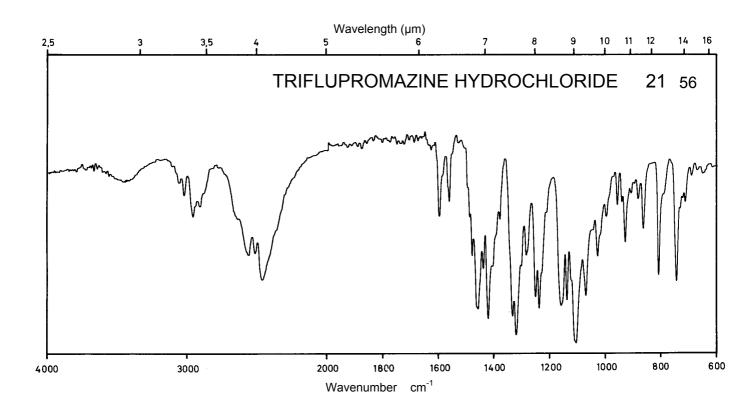



Name TRIFLUPERAZINE DIHYDROCHLORIDE

 M_{r} 480.4


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm 260 nm		306 nm 256 nm	310 nm 258 nm
E ^{1%} 1cm	77 697		73 628	69 593
3	3400 33460		3510 30180	3330 28500

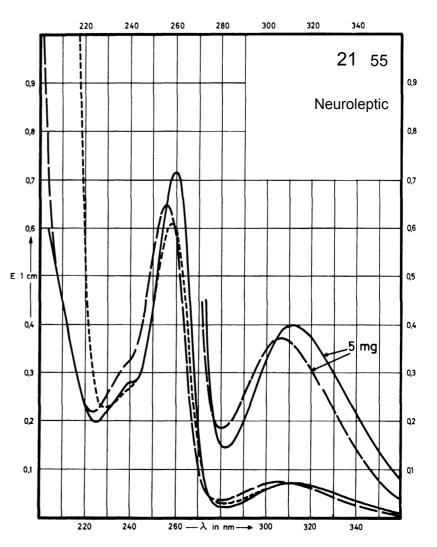

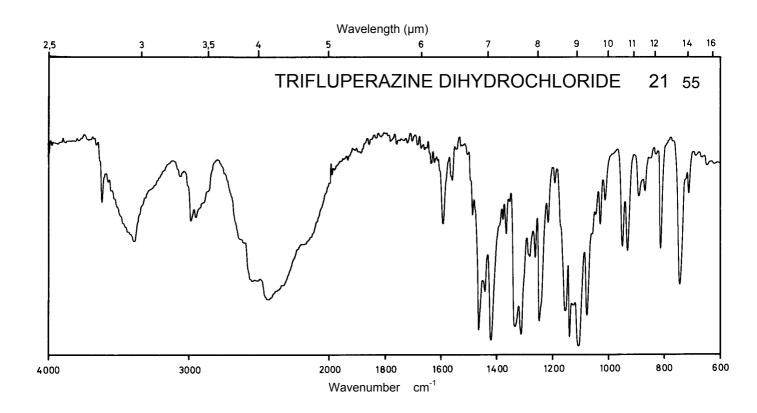

Name TRAZODONE HYDROCHLORIDE

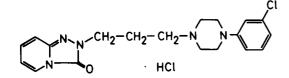



 M_r

408.3

Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	327 nm 252 nm		314 nm 247 nm	315 nm 247 nm
E ^{1%} 1cm	75 370		93 292	91 282
ε	3070 15100		3780 11920	3720 11530

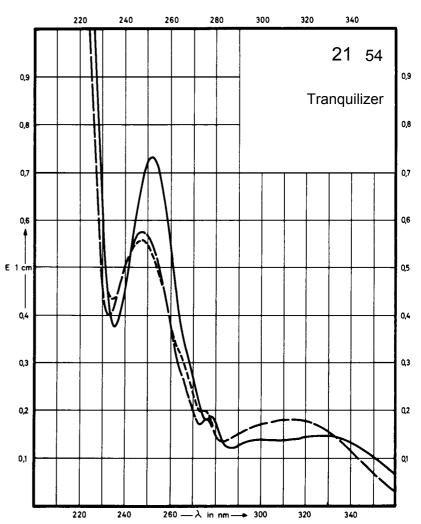

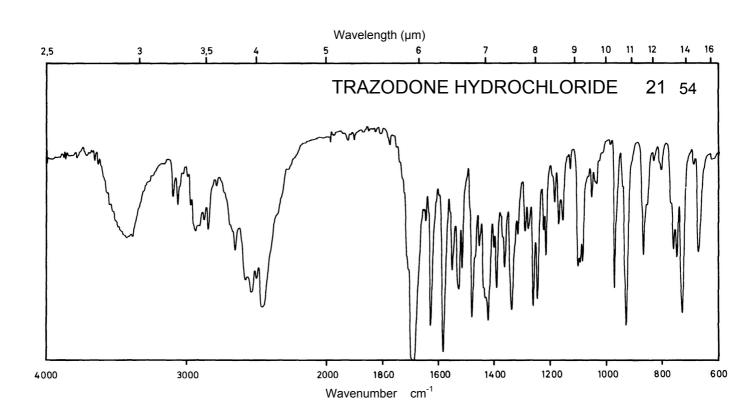



Name TIOTIXENE

 M_r 443.6


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm 260 nm 230 nm		308 nm 229 nm	310 nm 255 nm
E ^{1%} 1cm	180 346 830		182 816	170 342
3	8000 15360 36780		8060 36210	7530 15180

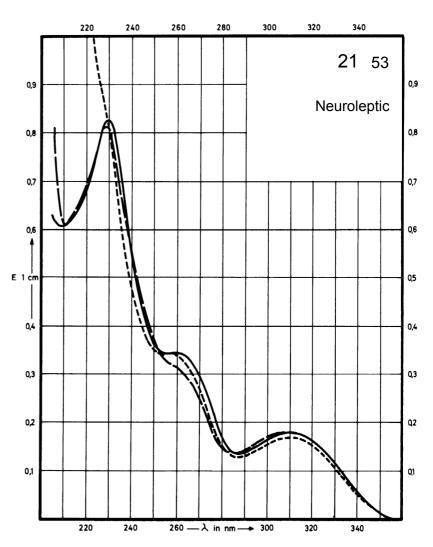

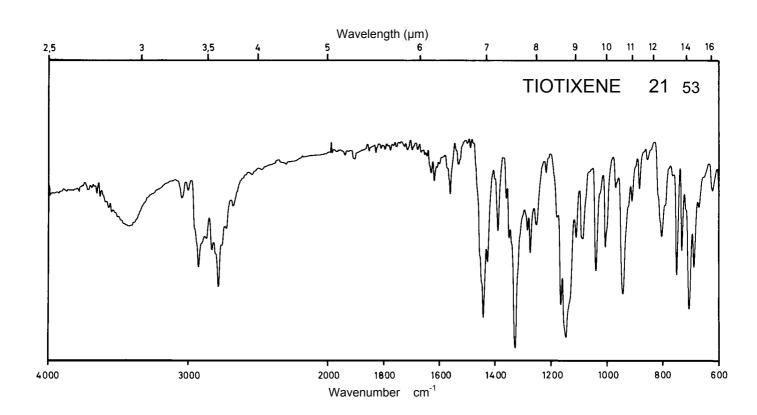



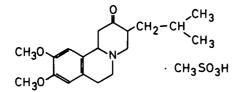



 M_r

413.5

Concentration 8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	284 nm		282 nm 228 nm	283 nm
E ^{1%} 1cm	89		85 180	91
ε	3670		3510 7440	3760

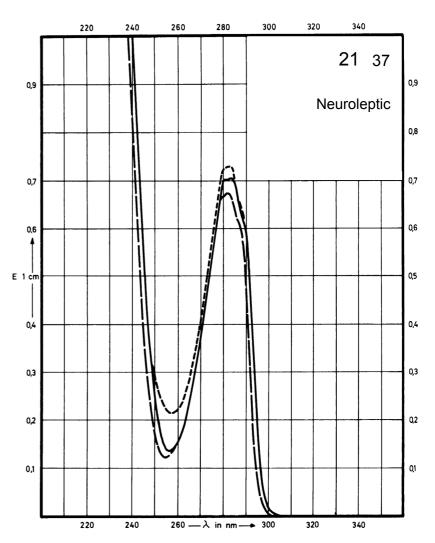

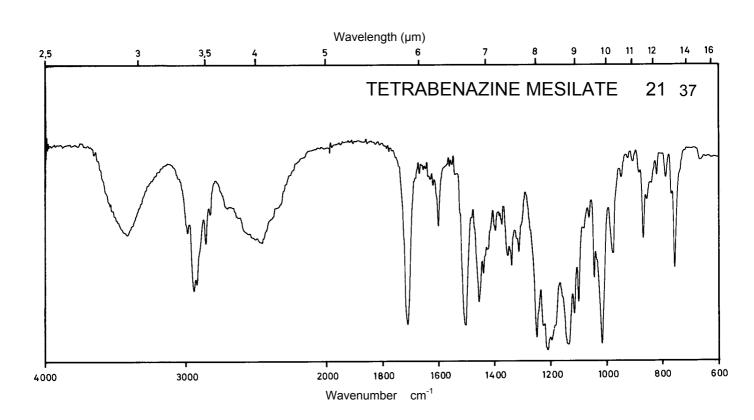


Name TRIMIPRAMINE MALEATE


 M_{r}

410.5


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	249 nm		250 nm	
E ^{1%} 1cm	238		211	
3	9780		8680	

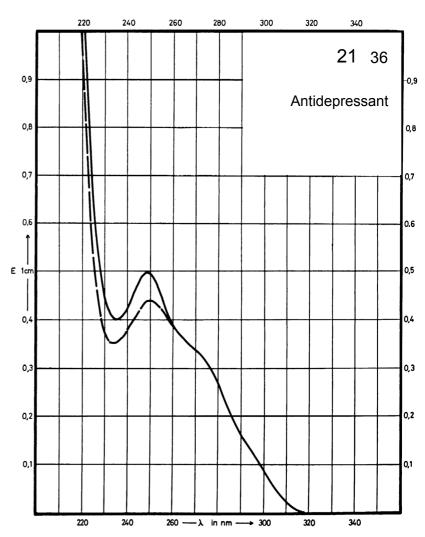

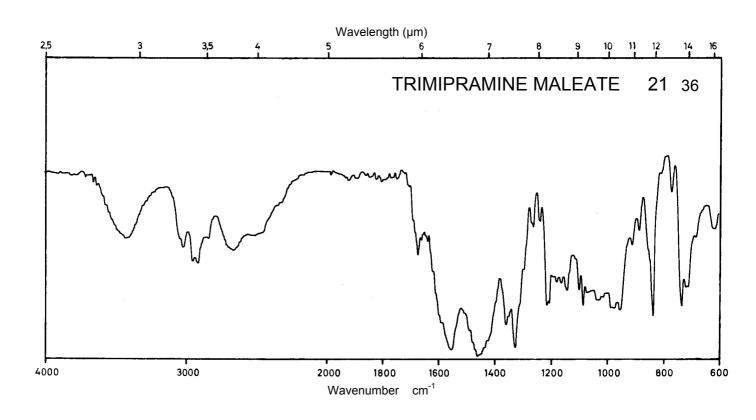



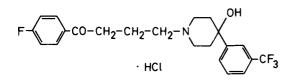



445.9

Concentration 2

2.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	244 nm		248 nm	
E ^{1%} 1cm	290		287	
ε	12950		12770	

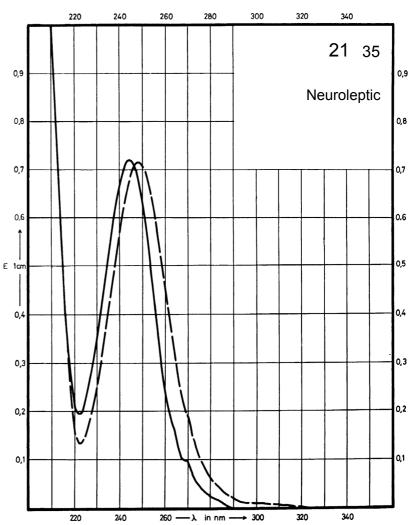

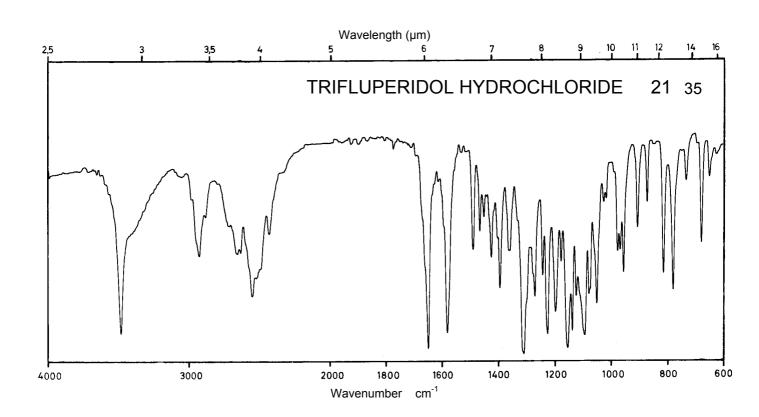

Name TRANYLCYPROMINE SULFATE


$$\begin{bmatrix} \mathsf{NH}_3 \\ \mathsf{CH} \\ \mathsf{CH} - \mathsf{CH}_2 \end{bmatrix}_2^* \quad \mathsf{SO}_4^{--}$$


 M_{r}

364.5

Concentration 30 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm 265 nm		271 nm 264 nm	274 nm 267 nm
E 1% 1cm	12.1 16.6		11.5 16.0	20.8 27.8
ε	443 603		419 585	758 1010

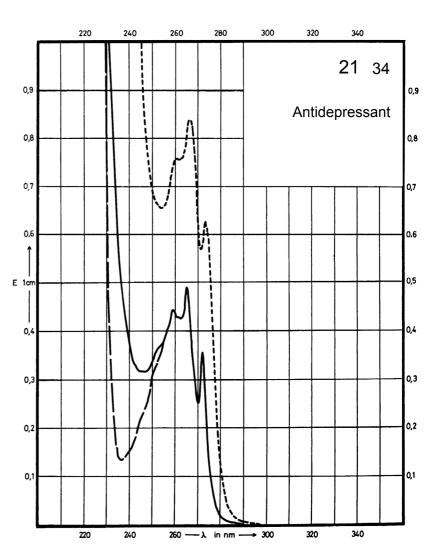

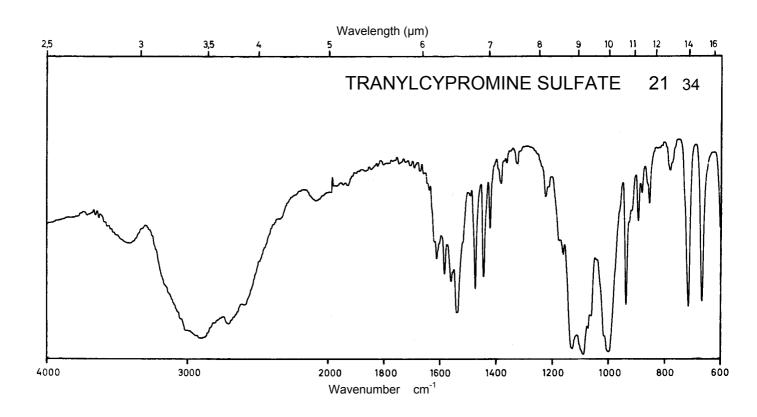


Name THIORIDAZINE HYDROCHLORIDE


 M_{r}

407.0


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm 263 nm		310 nm 262 nm	
E ^{1%} 1cm	112 936		107 932	
ε	4560 38010		4350 37850	

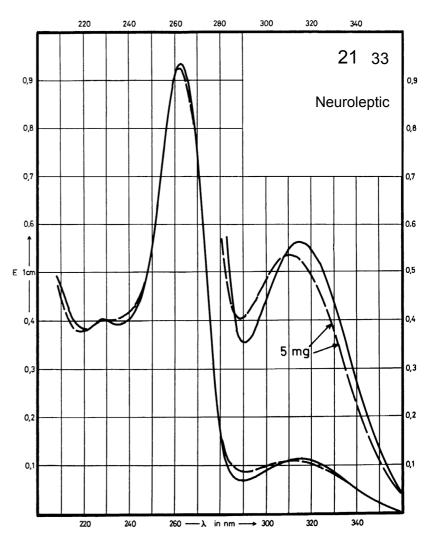

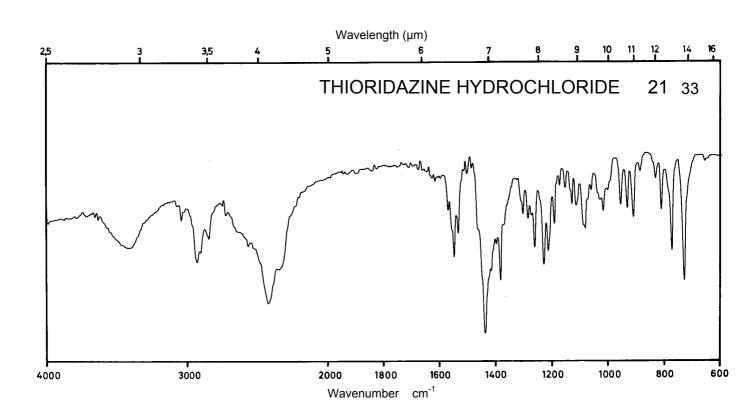




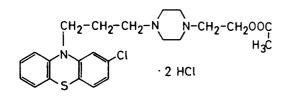


518.9


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	308 nm 256 nm		305 nm 255 nm	308 nm 257 nm
E ^{1%} 1cm	81 675		78 640	70 542
3	4200 35030		4050 33210	3630 28120

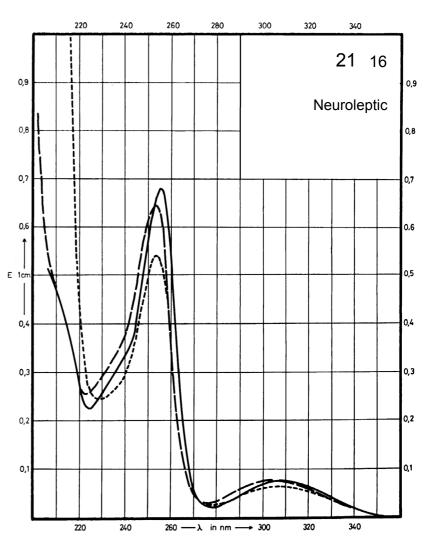

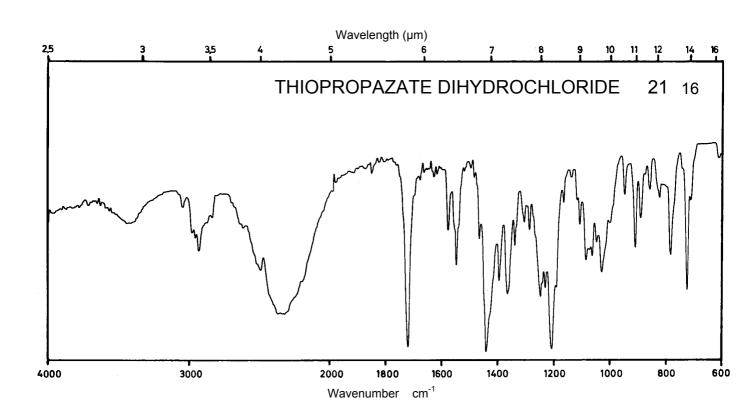



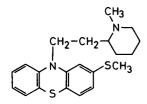



370.6

Concentration 0.5 mg / 100 ml 5 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH*
Maximum of absorption	315 nm 263 nm		310 nm 262 nm	315 nm 264 nm
E ^{1%} 1cm	129 1072		118 1028	129 1095
ε	4770 39720		4360 38110	4770 40620

* 1 M NaOH + Methanol (1+9)

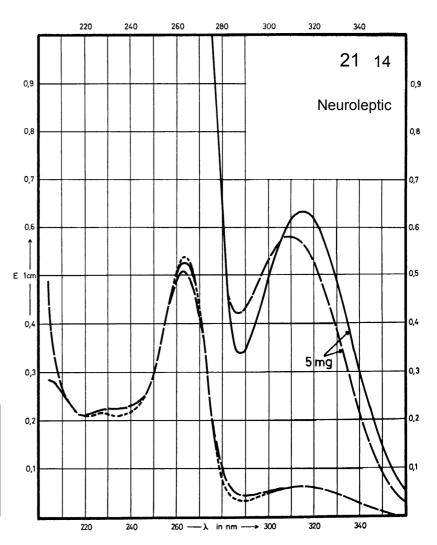

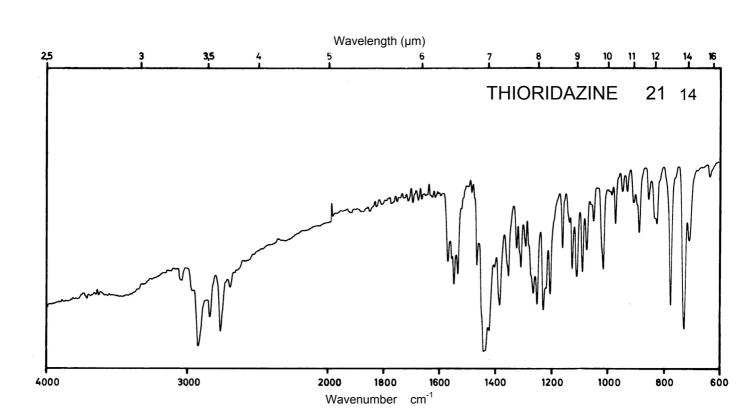




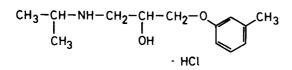


 M_r 259.8


Concentration 2.5 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	279 nm 272 nm 218 nm		278 nm 272 nm	278 nm 272 nm
E ^{1%} 1cm	52 55 281		45 51	46 52
ε	1360 1430 7340		1180 1320	1200 1360

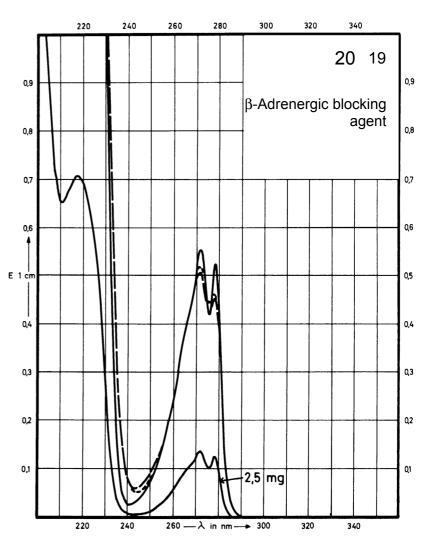

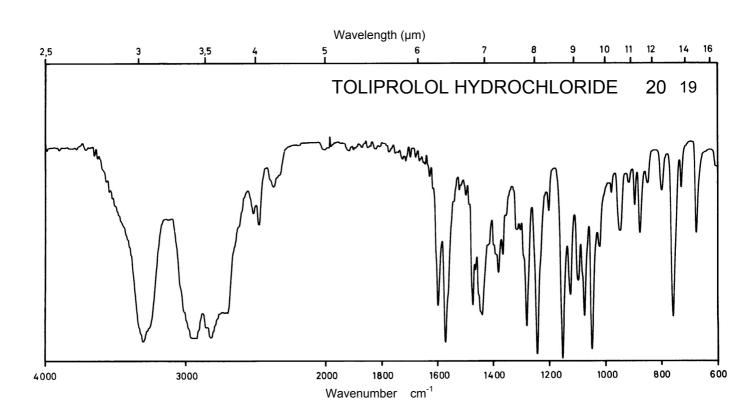



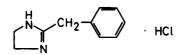



 M_r

196.7


Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm 258 nm	264 nm 257 nm		
E ^{1%} 1cm	7.5 11.4	7.5 11.3		
ε	148 225	148 222		

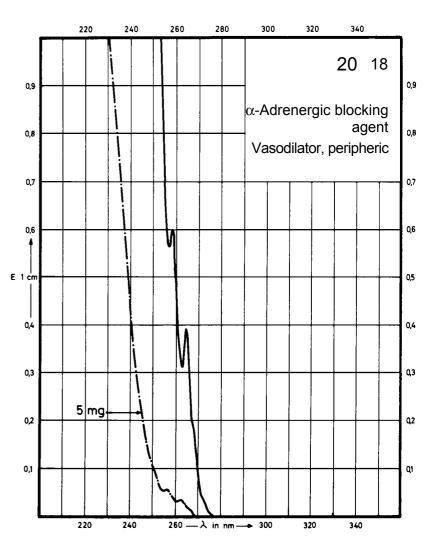

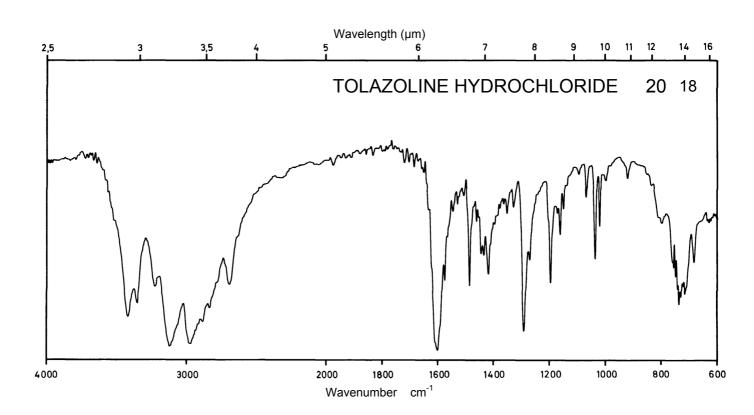




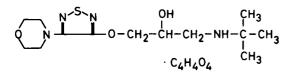


 M_r 432.5


Concentration 3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	299 nm		295 nm	296 nm
E ^{1%} 1cm	210		197	206
3	9070		8540	8920

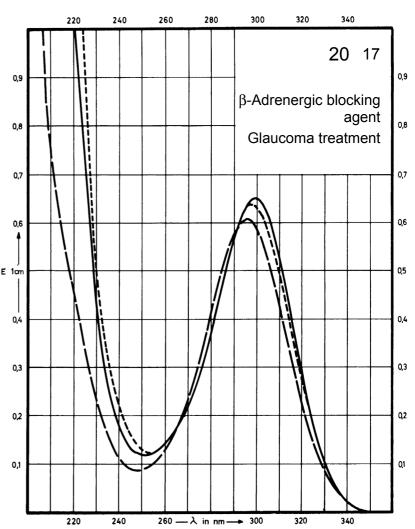

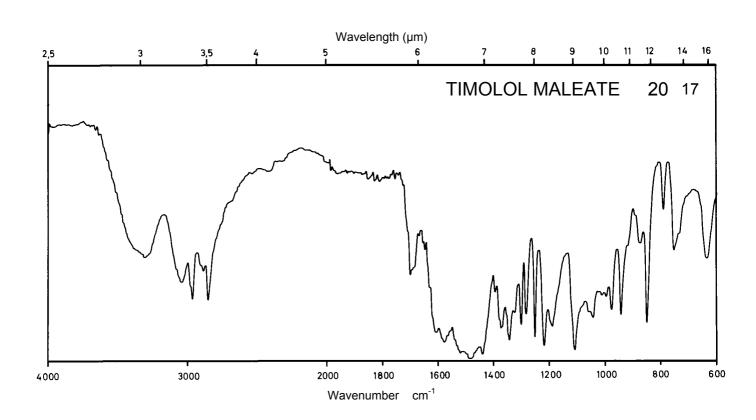



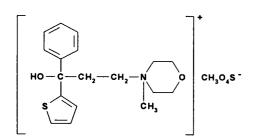



 M_r

429.6

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	233 nm	233 nm	233 nm	
E ^{1%} 1cm	220	208	210	
ε	9400	8900	9000	

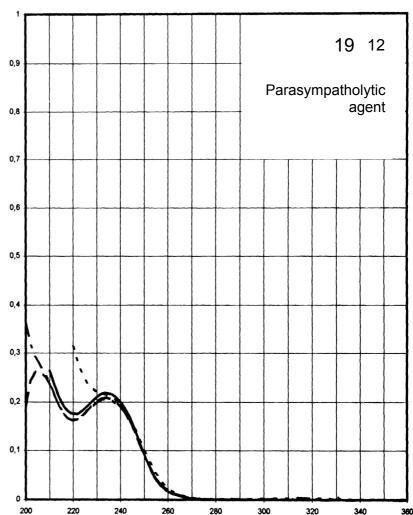

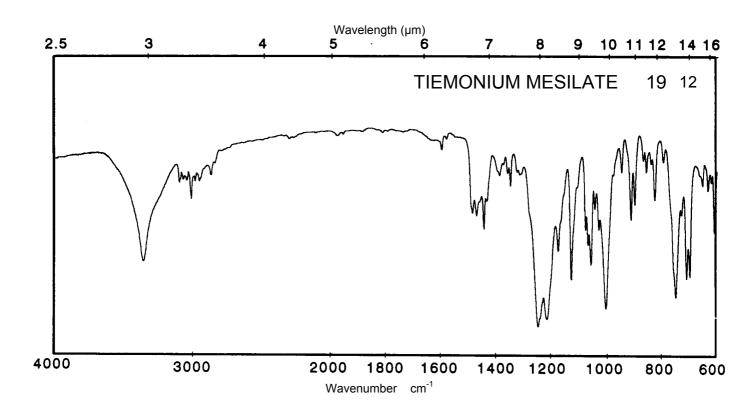

Name TIEMONIUM IODIDE


 M_{r}


445.4

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	221 nm	227 nm	227 nm	
E ^{1%} 1cm	492	488	486	
3	21900	21800	21700	

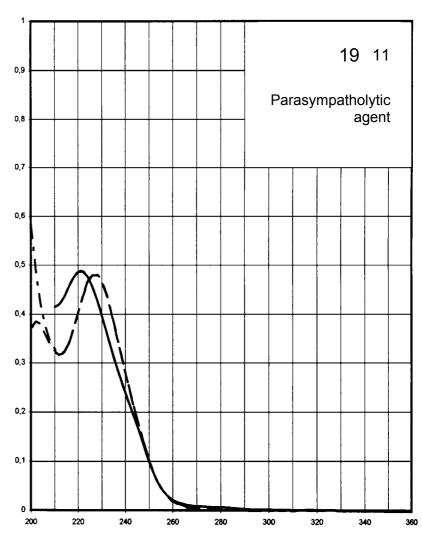

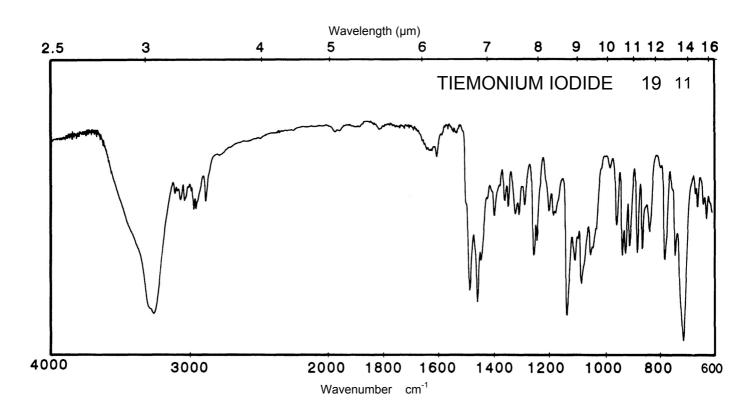




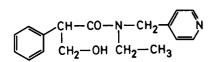


284.4

Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	256 nm		253 nm	256 nm
E ^{1%} 1cm	80		174	84
ε	2280		5000	2400

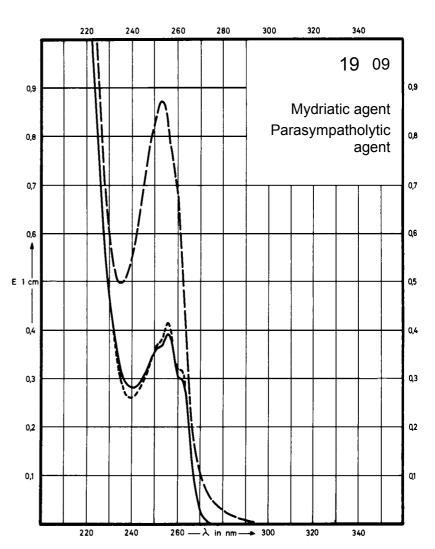

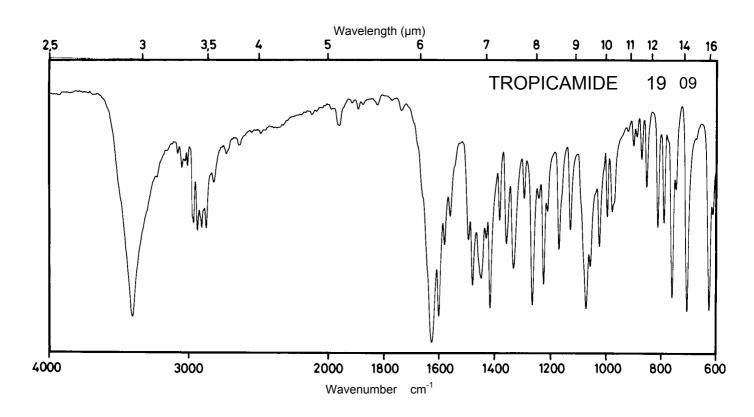


Name TOLYCAINE HYDROCHLORIDE


 M_r

314.8


Concentration 2 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	286 nm	284 nm	284 nm	Decom- position observed
E ^{1%} _{1cm}	61	58	58	
ε	1910	1830	1830	

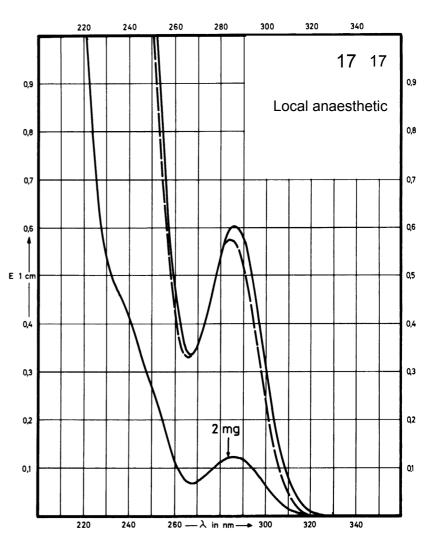

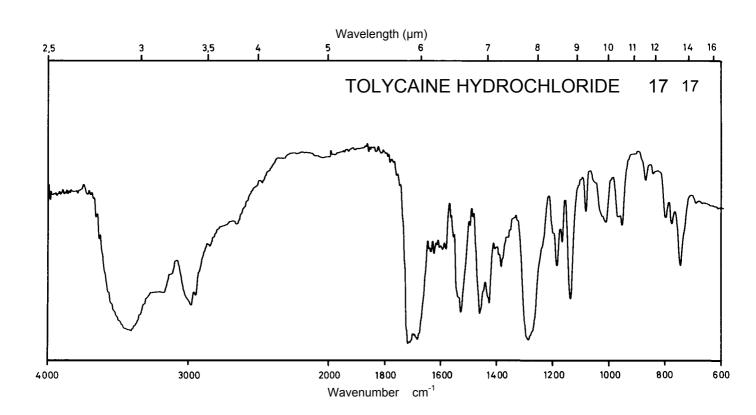



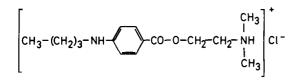



 M_r

300.8

Concentration 1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	310 nm 226 nm		310 nm 278 nm 227 nm	303 nm
E ^{1%} 1cm	967 242		49 46 472	666
ε	29090 7280		1470 1380 14200	20030

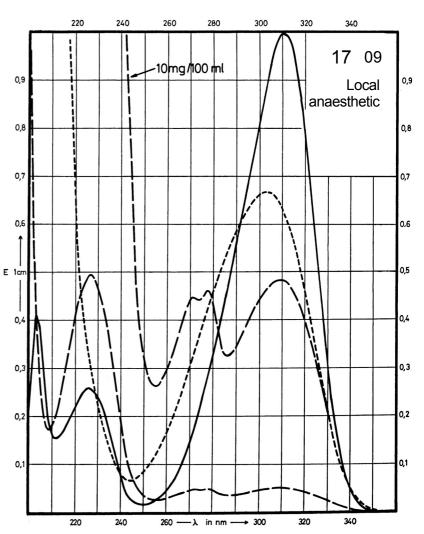

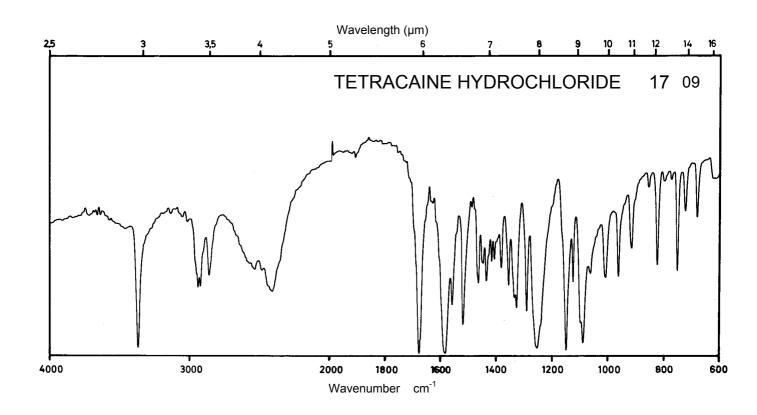

Name THIOBUTABARBITAL SODIUM


$$\begin{bmatrix} 0 & H & S \\ H_5C_2 & HC & N \\ H_3C & 0 \end{bmatrix}$$
 Na⁺


 M_{r}

250.3

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			287 nm 238 nm	304 nm
E ^{1%} 1cm			711 288	809
ε			17800 6960	20240

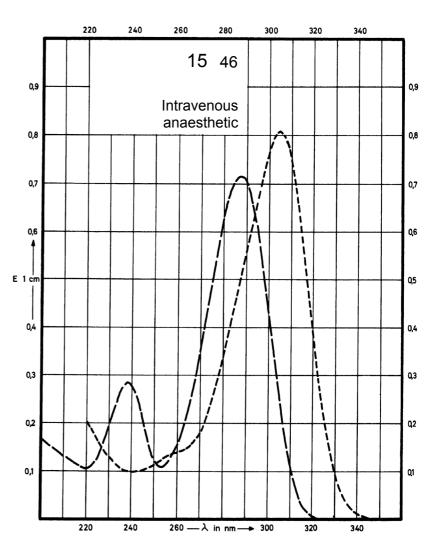

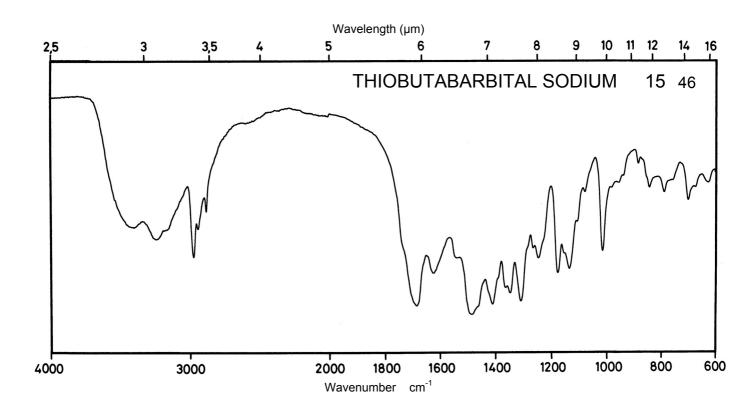





H₃C O CH₃

 M_r 143.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	227 nm			
E ^{1%} 1cm	9.0			
ε	130			

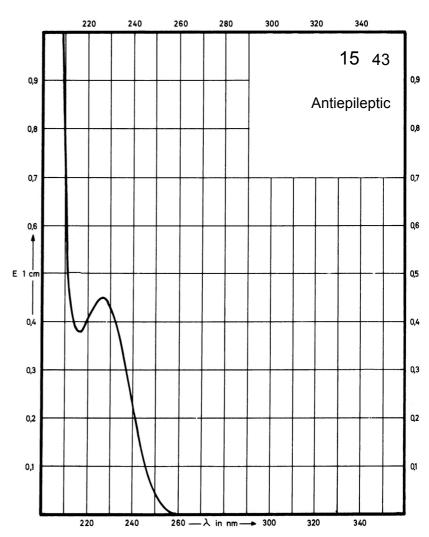

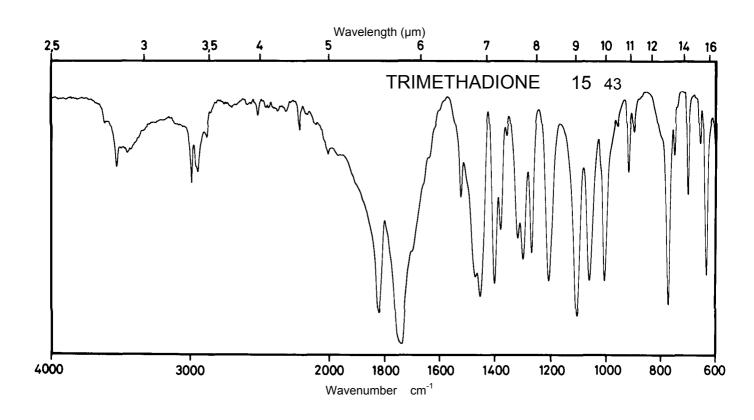



Name TRIAZOLAM

CH₃ N N


 M_r 343.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	221 nm	221 nm	Decom- position observed	
E ^{1%} 1cm	1200	1190		
ε	41180	40840		

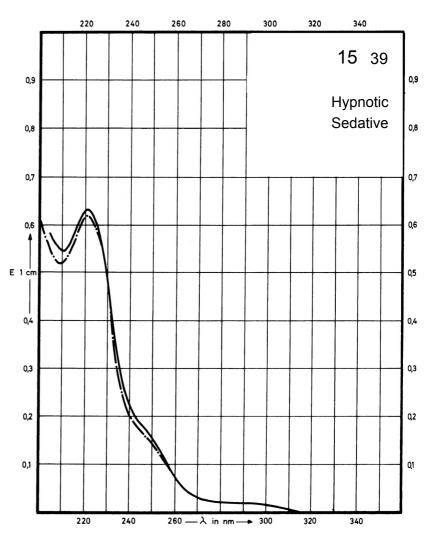

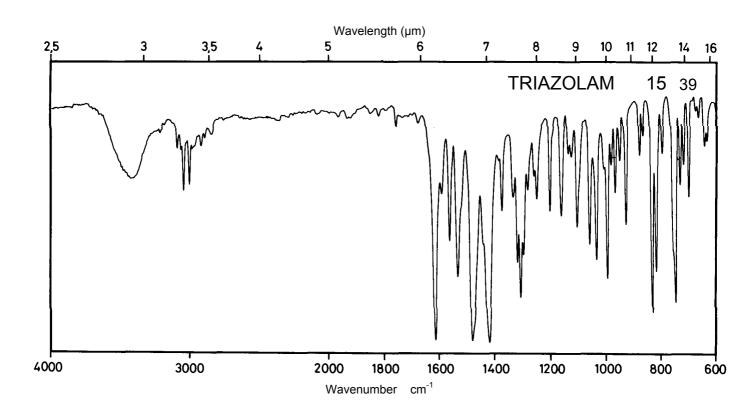


Name THALIDOMIDE

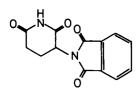


 M_r 258.2


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	294 nm 240 nm		300 nm 241 nm	
E ^{1%} 1cm	74 358		94 385	
3	1920 9240		2430 9950	

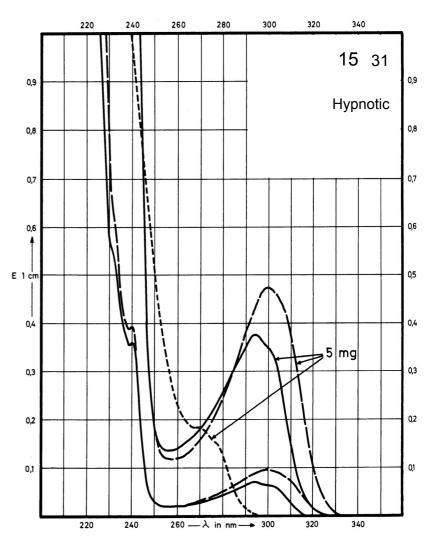

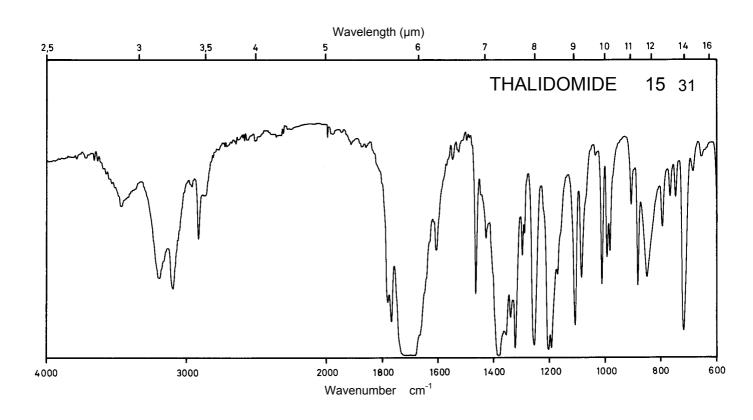


Name THIOPENTAL SODIUM

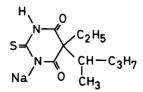


 M_{r}


264.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm 237 nm		287 nm 238 nm	305 nm
E ^{1%} 1cm	792 320		817 320	933
ε	20930 8460		21600 8460	24660

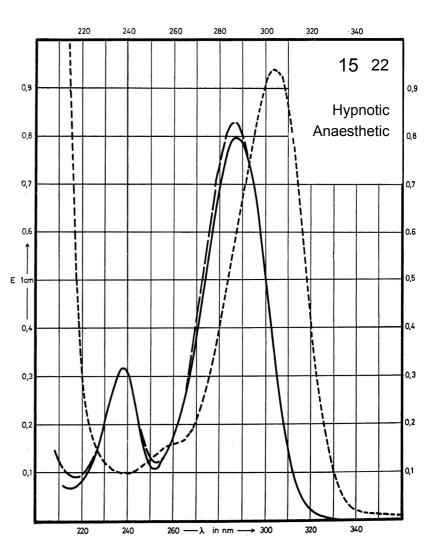

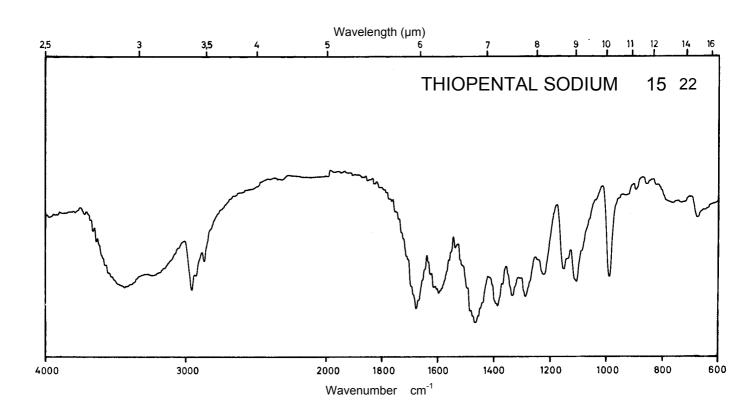





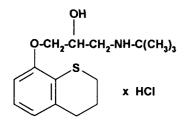

 M_{r} 331.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	293 nm 257 nm 218 nm	289 nm 254 nm 216 nm	289 nm 253 nm 216 nm	291 nm 254 nm
E ^{1%} 1cm	64 286 817	59 258 815	59 258 816	59 240
3	2100 9500 27100	1960 8500 27000	1970 8600 27100	1970 8000

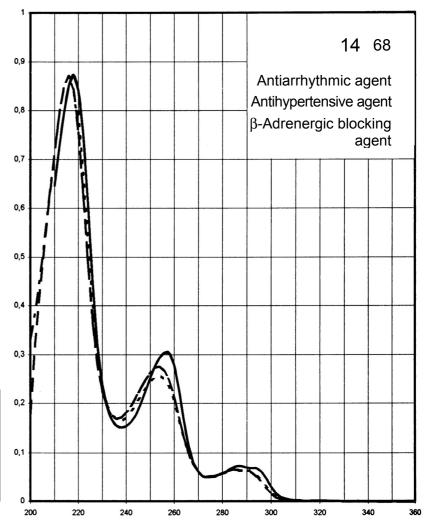

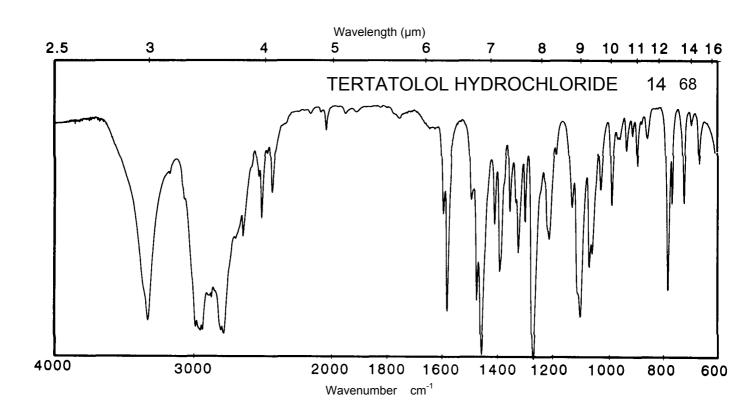


Name TRIMETAZIDINE DIHYDROCHLORIDE


 M_{r}

339.3

Concentration 2.5 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm		268 nm	270 nm
E ^{1%} 1cm	23.2		22.4	24.3
3	790		760	830

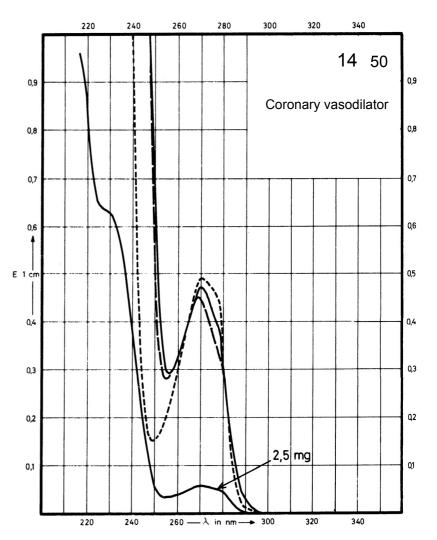

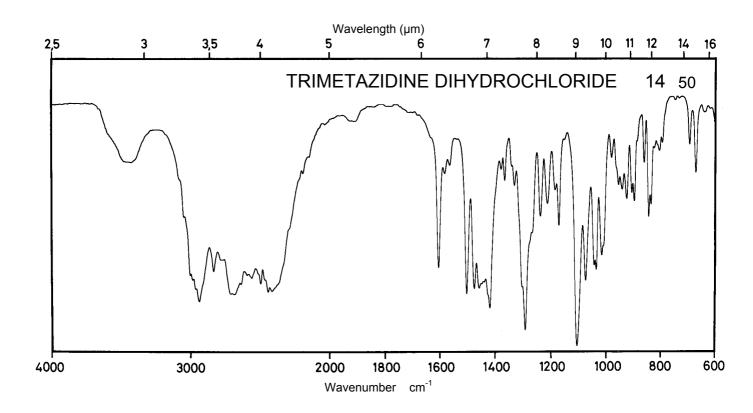



Name TOCAINIDE HYDROCHLORIDE

 M_{r}

228.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm 264 nm	272 nm 263 nm	272 nm 263 nm	
E ^{1%} 1cm	12.0 15.5	15.2 18.5	15.2 18.5	
ε	275 355	348 425	348 425	

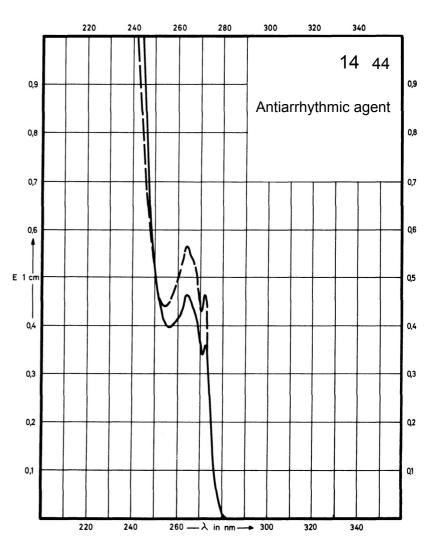

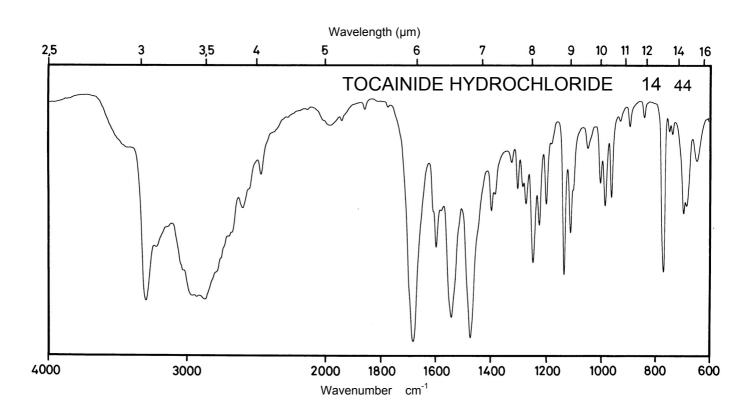



Name THEOPHYLLINE

H₃C N H₂O . H₂O


 M_r 198.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm		270 nm	275 nm
E ^{1%} 1cm	501		538	632
ε	9930		10660	12530

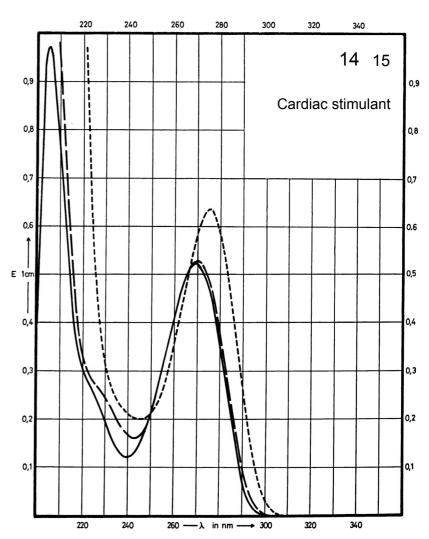

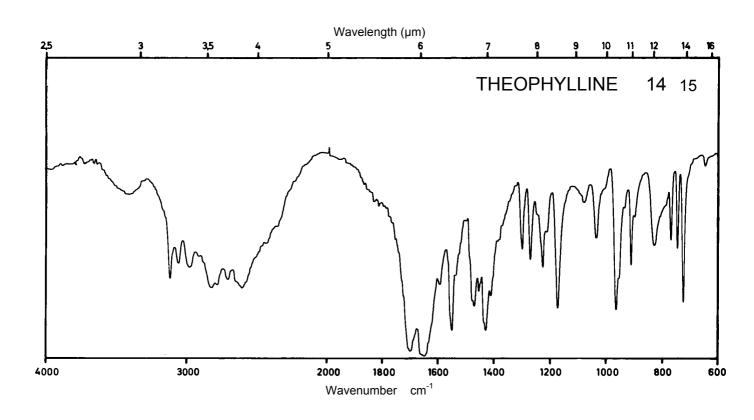





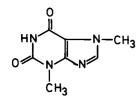


180.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		271 nm	273 nm
E ^{1%} 1cm	532		561	574
ε	9590		10110	10340

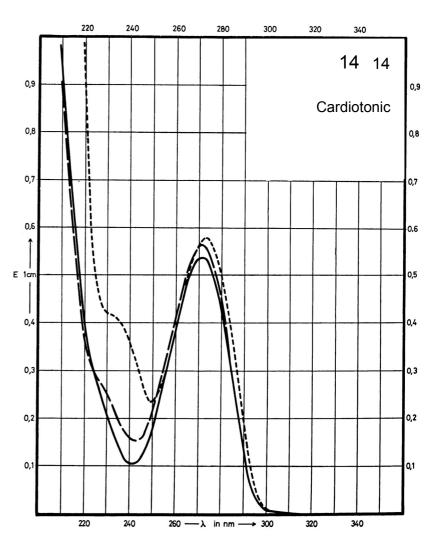

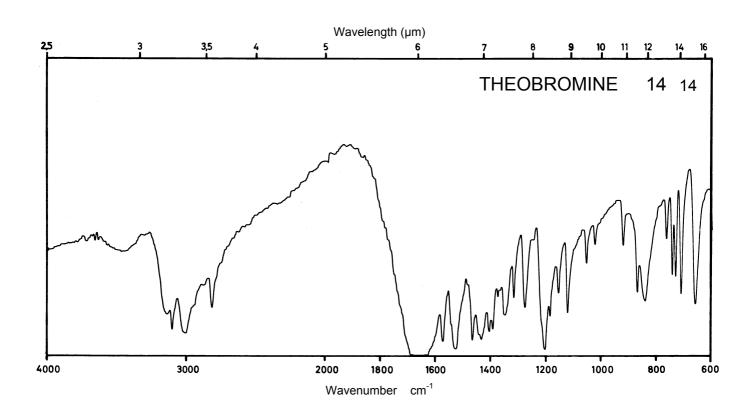





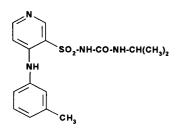


 M_r 348.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm	290 nm	286 nm	284 nm
E ^{1%} 1cm	344	378	347	364
ε	12000	13200	12100	12700

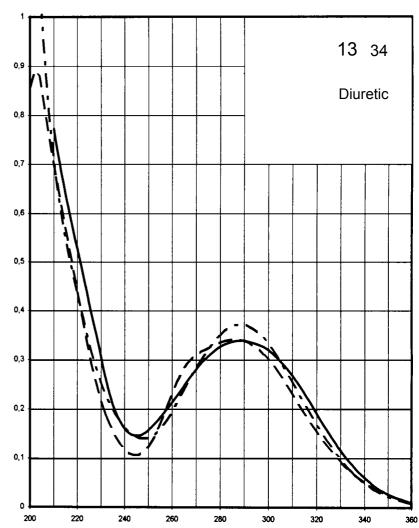

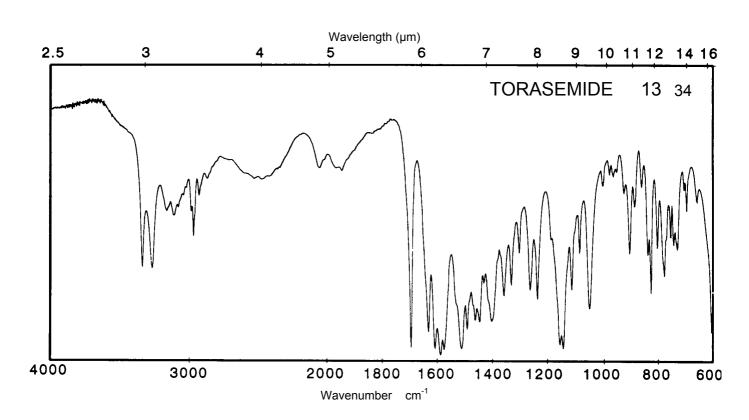


Name TRICHLORME-THIAZIDE

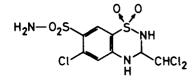


 M_r 380.7

Concentration 0.5 mg / 100 ml 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	312 nm 267 nm 225 nm	312 nm 268 nm 225 nm	312 nm 268 nm 225 nm	317 nm 269 nm
E 1% 1cm	78 560 1120	82 537 1120	82 546 1140	70 438
ε	3000 21300 42600	3100 20500 42600	3100 20800 43400	2660 16700

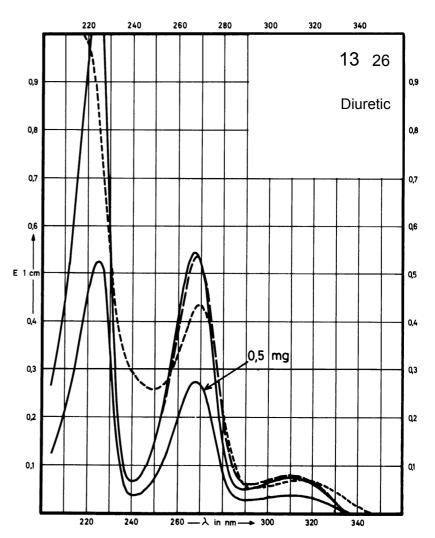

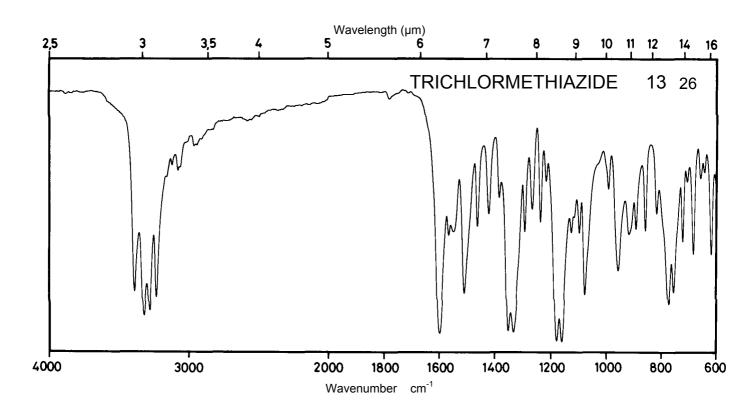



Name TIENILIC ACID

 M_{r}

331.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	295 nm 268 nm		304 nm 274 nm	305 nm 274 nm
E ^{1%} 1cm	387 345		388 342	397 339
ε	12810 11420		12840 11330	13140 11230

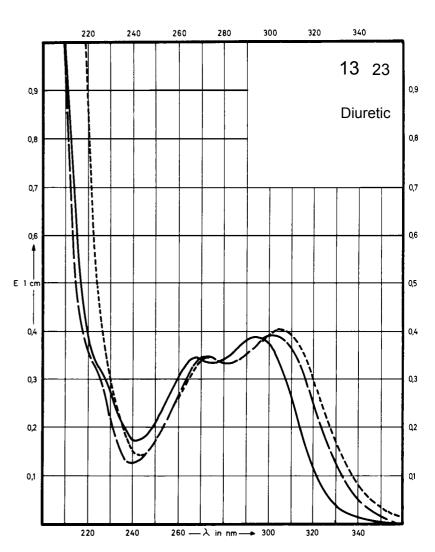

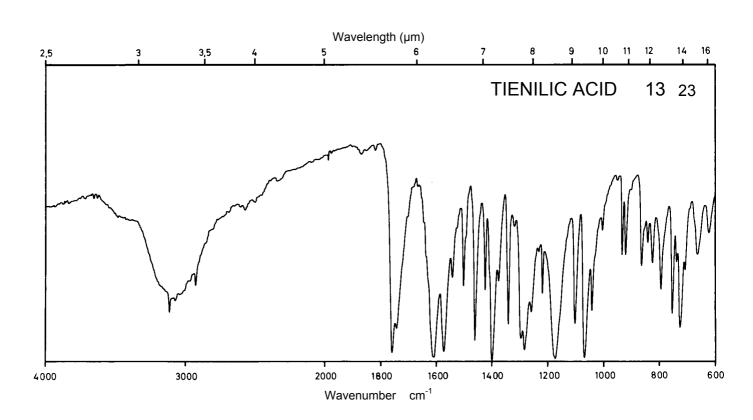



Name TIZOLEMIDE HYDROCHLORIDE

 M_{r}


372.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	282 nm	282 nm	282 nm	282 nm
	274 nm	274 nm	274 nm	273 nm
E ^{1%}	27.4	34.6	34.6	20.4
1cm	29.8	37.7	37.5	27.4
ε	1020	1290	1290	760
	1110	1400	1400	1020

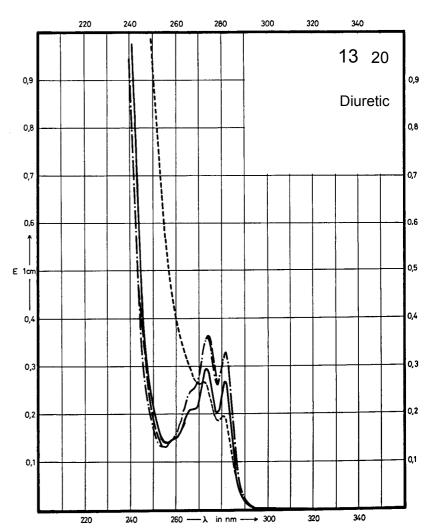

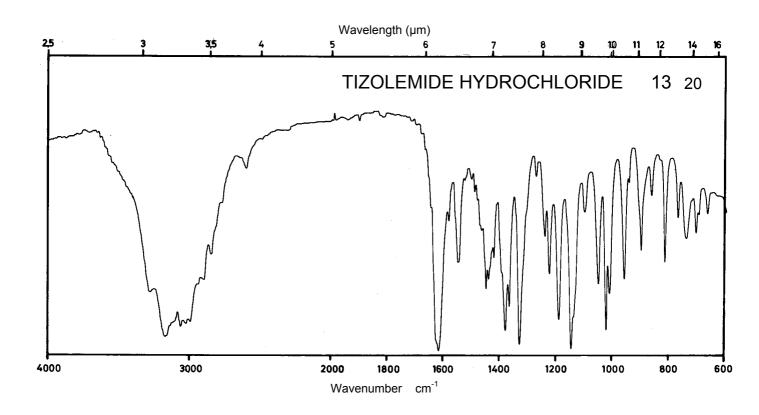



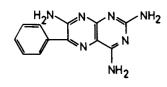



253.3

Concentration


0.9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	266 nm		357 nm 250 nm	270 nm
E ^{1%} 1cm	568		855 634	559
ε	14390		21660 16060	14160

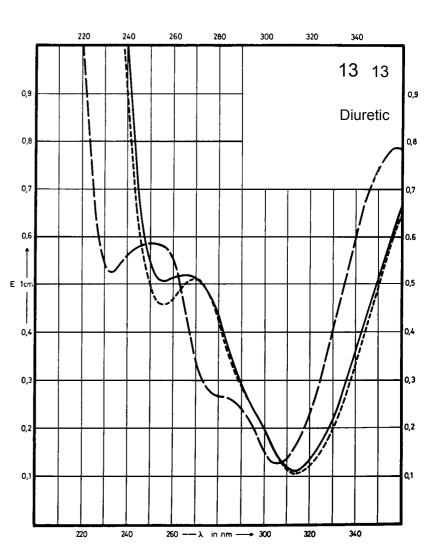

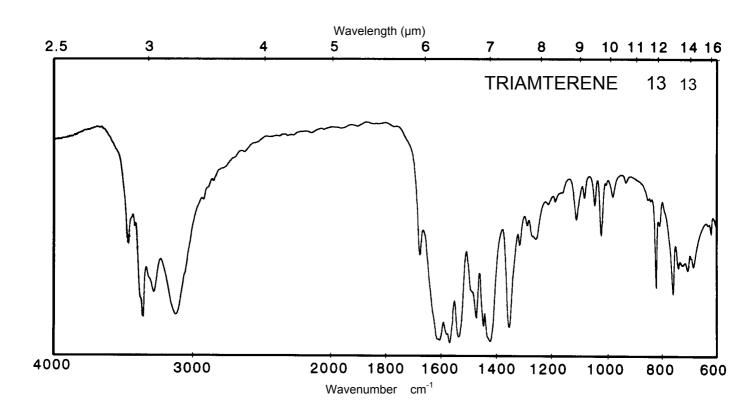





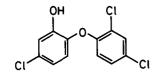

289.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	281 nm		279 nm	292 nm
E ^{1%} 1cm	180		140	230
3	5200		4060	6650

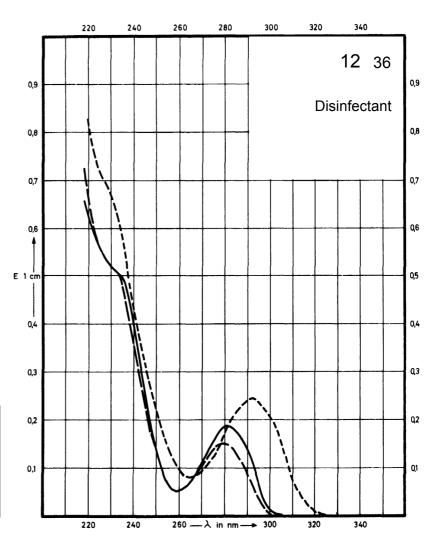

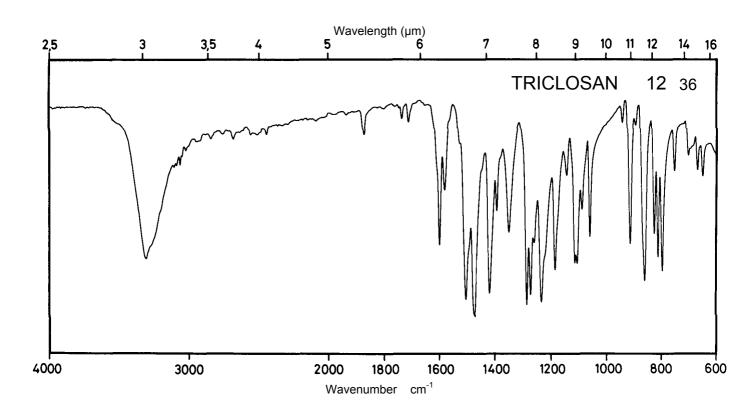


Name TOSYLCHLORAMIDE SODIUM


 M_{r}

281.7

Concentration 2 mg / 100 ml 21 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	220 nm		274 nm 228 nm	
E ^{1%} 1cm	352		28 473	
ε	9920		790 13320	

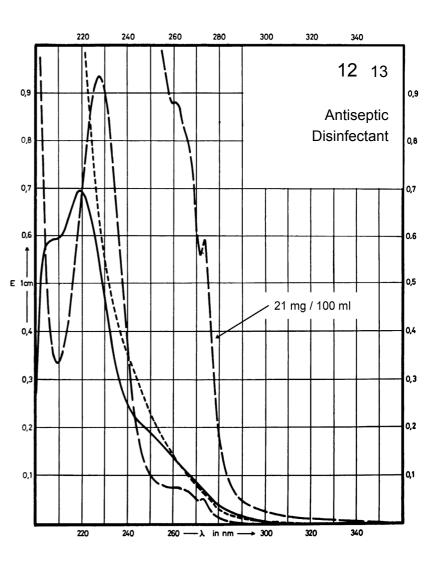

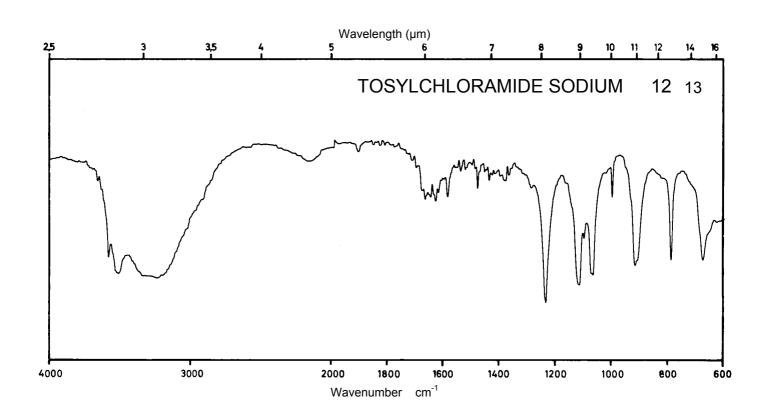



Name THYMOL

 H_3C CH CH_3


 M_r 150.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	276 nm		272 nm	292 nm
E ^{1%} 1cm	158		137	257
ε	2370		2060	3860

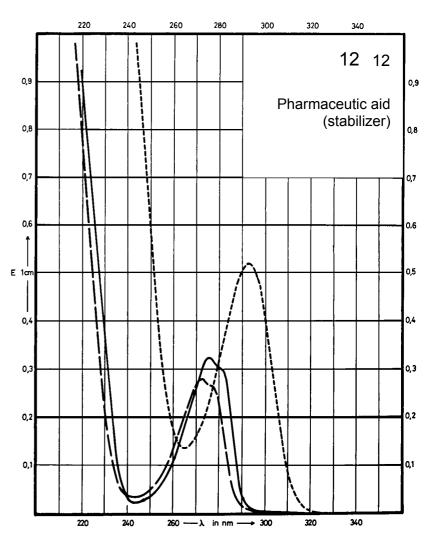

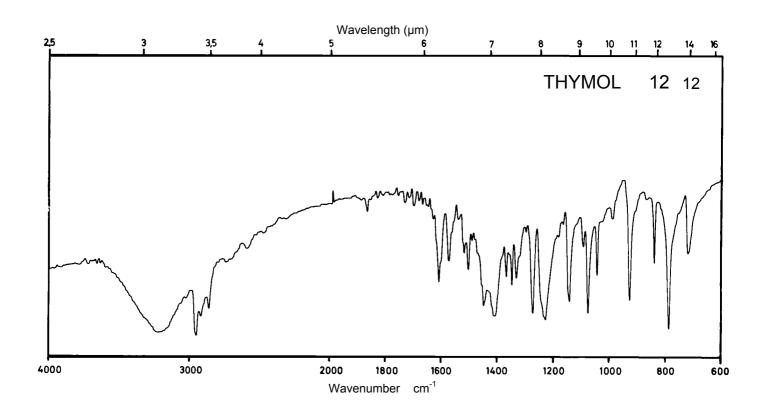




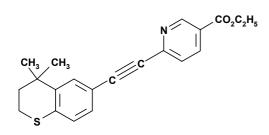


 M_{r}


351.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	351 nm 259 nm			
E ^{1%} 1cm	988 437			
3	34700 15400			

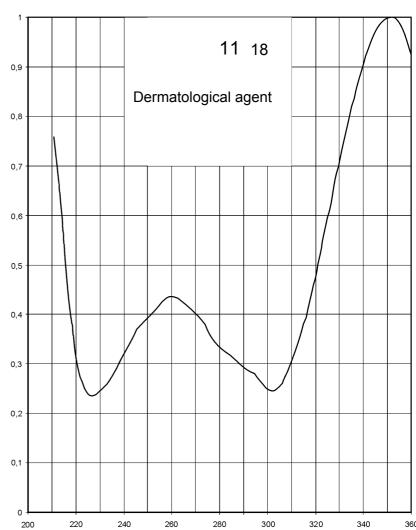

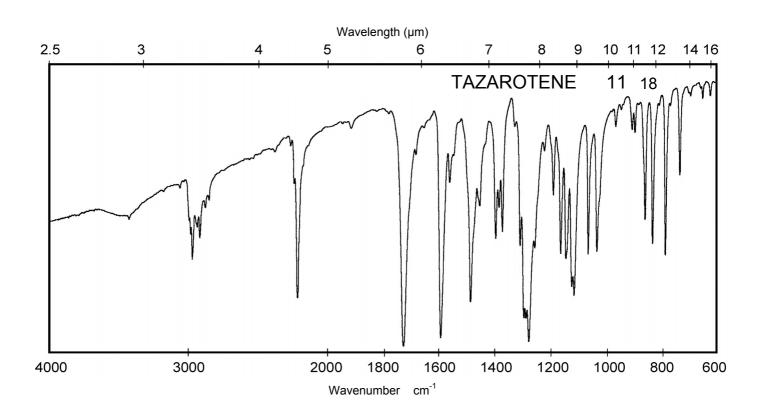





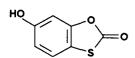



168.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 233 nm	285 nm 231 nm	285 nm 231 nm	Decom- position observed
E ^{1%} 1cm	255 587	260 551	262 552	
ε	4300 9900	4600 9300	4400 9300	

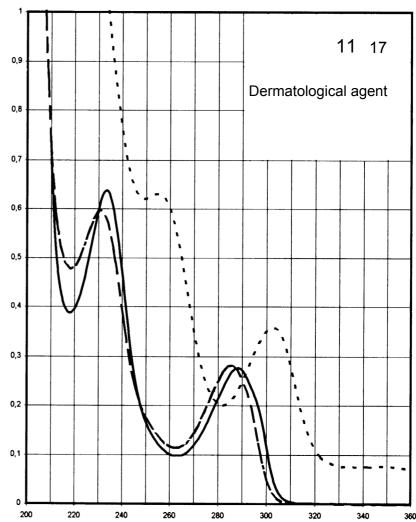

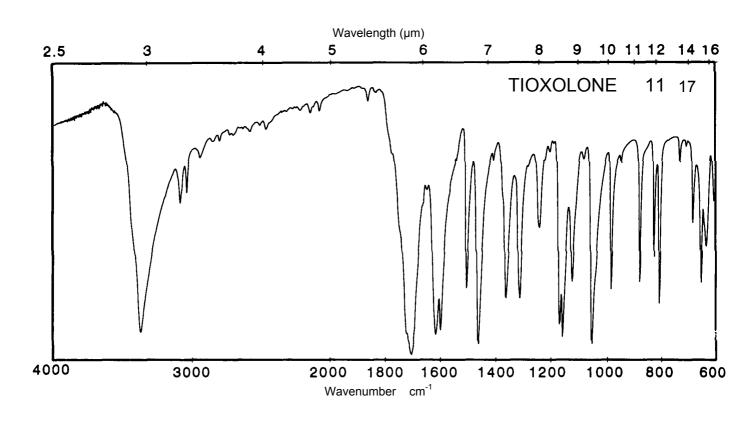

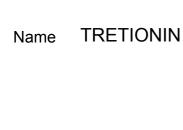

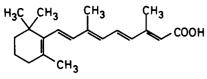



300.4

Concentration

0.32 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	350 nm			
E ^{1%} 1cm	1480			
ε	44500			

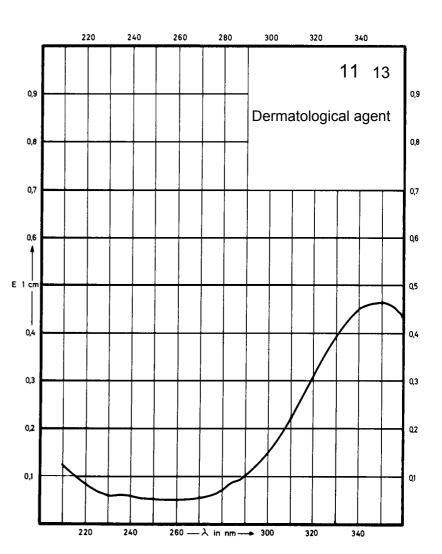

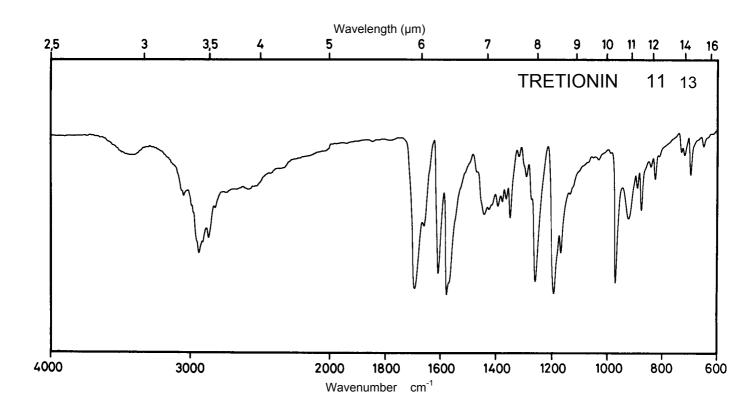



Name TRIAMCINOLONE HEXACETONIDE

 M_{r}


532.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm			
E ^{1%} 1cm	291			
ε	15480			

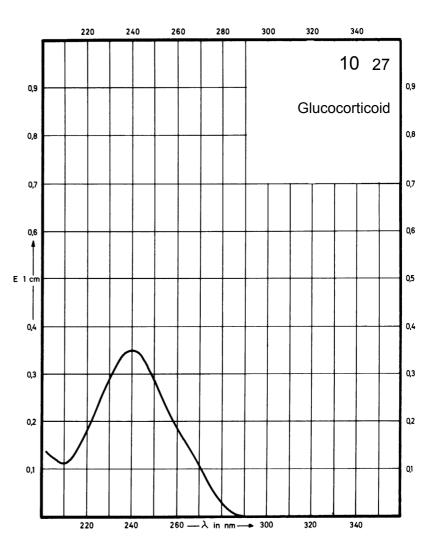

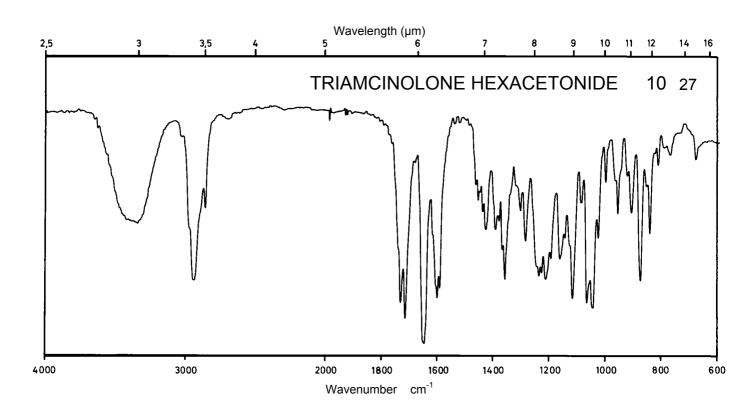


Name TRIAMCINOLONE 16α, 21-DIACETATE

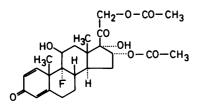


 $\,M_r\,$


478.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm	242 nm		
E ^{1%} 1cm	312	315		
ç	14910	15080		

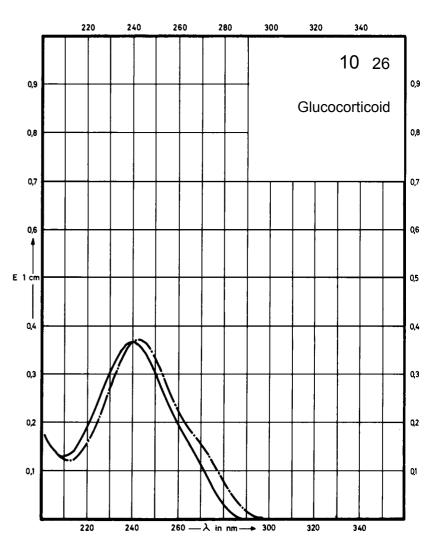

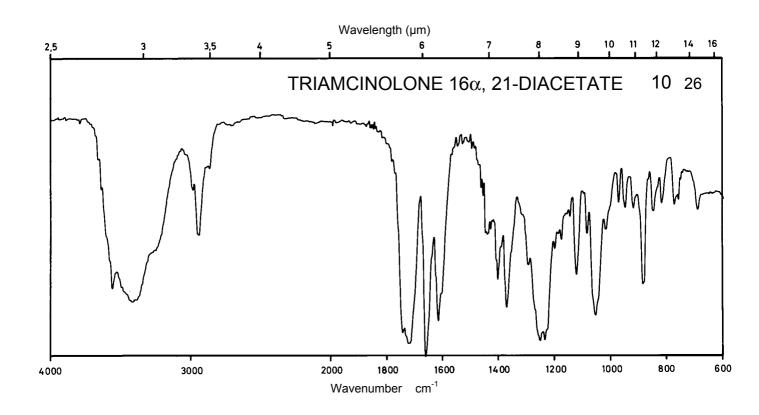



Name TRIAMCINOLONE 16α , 17α -ACETONIDE

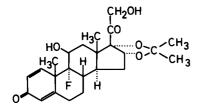


 M_r 434.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	239 nm	243 nm		
E ^{1%} 1cm	347	355		
3	15060	15430		

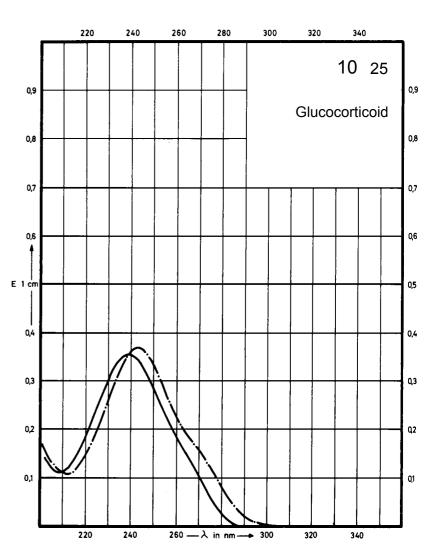

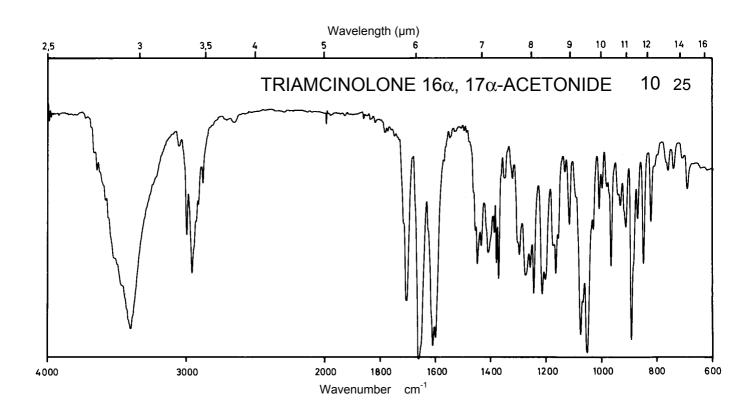





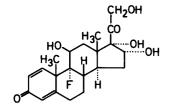


 M_r 394.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	240 nm	242 nm		
E ^{1%} 1cm	384	390		
ε	15140	15380		

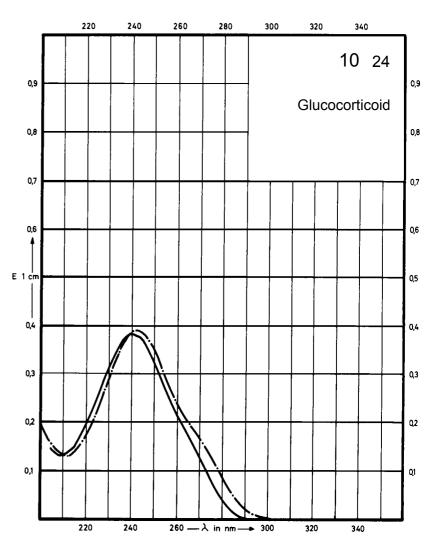

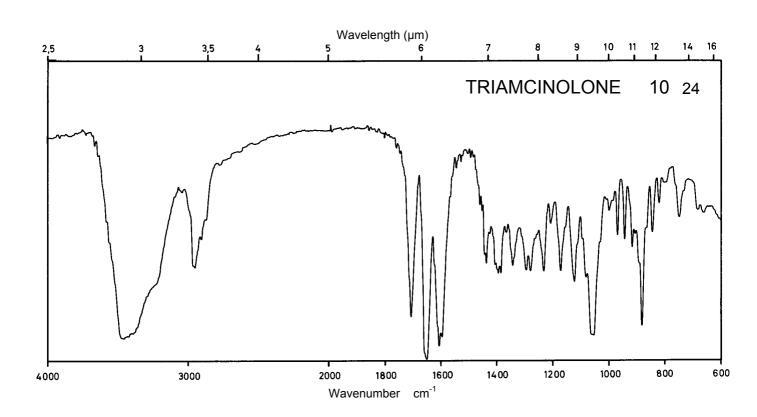


Name TRENBOLONE 17β -ACETATE

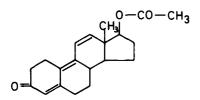


 M_{r}


312.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	339 nm 236 nm		351 nm 238 nm	352 nm 238 nm
E 1% 1cm	942 199		863 220	872 220
ε	29430 6220		26960 6870	27240 6870

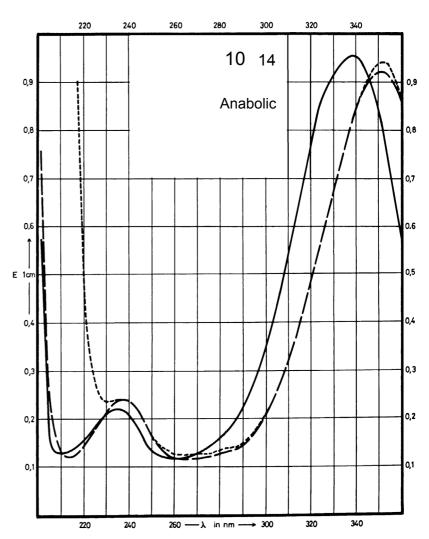

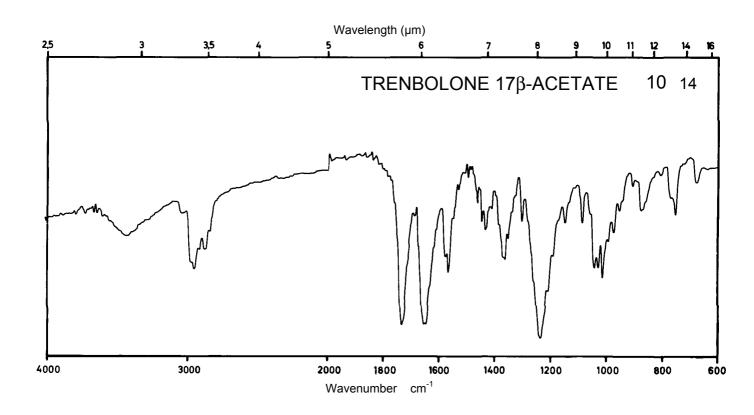





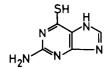


 M_{r} 167.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			347 nm 257 nm	320 nm 228 nm
E ^{1%} 1cm			1245 490	1015 1060
ε			20800 8200	17000 17700

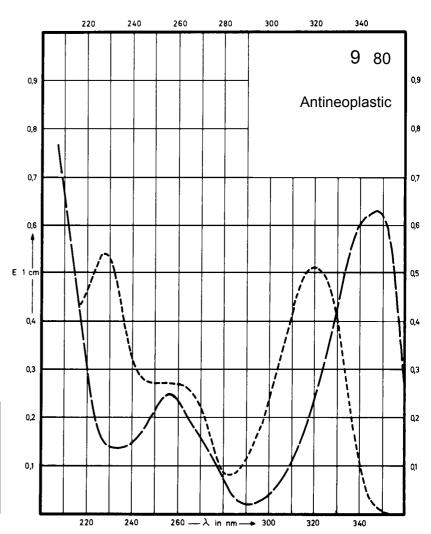

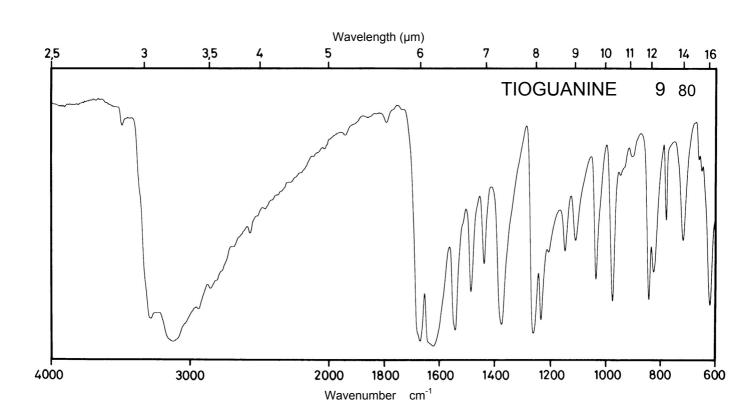


Name TINIDAZOLE

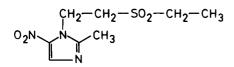


 $\,M_r\,$


247.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	311 nm 230 nm	318 nm 231 nm	277 nm	Decom- position observed
E ^{1%} 1cm	354 148	369 143	245	
3	8750 3630	9120 3530	6070	

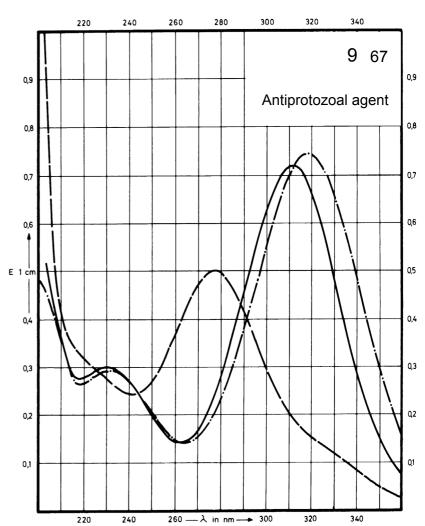

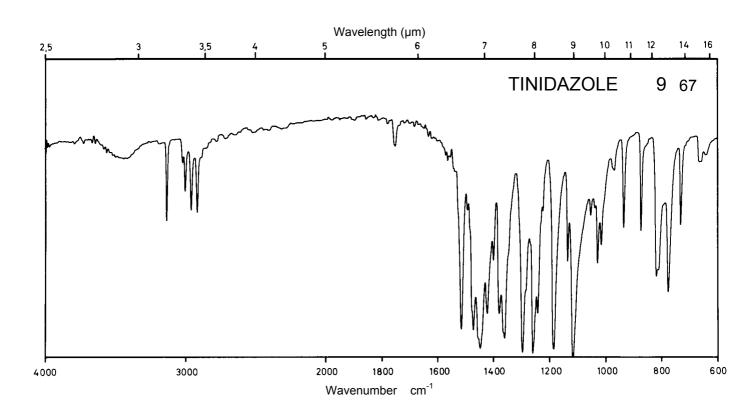





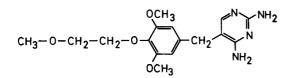


 M_r 334.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm	279 nm	270 nm	288 nm
E ^{1%} 1cm	202	171	178	219
3	6770	5730	5960	7320

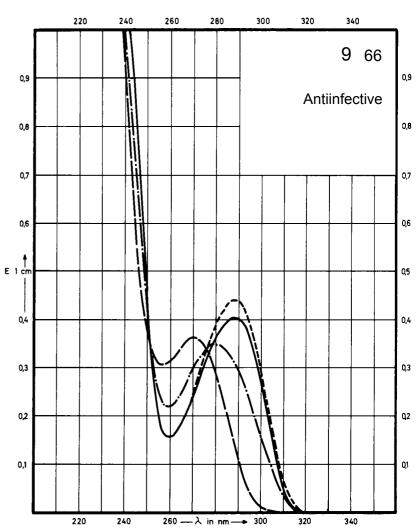

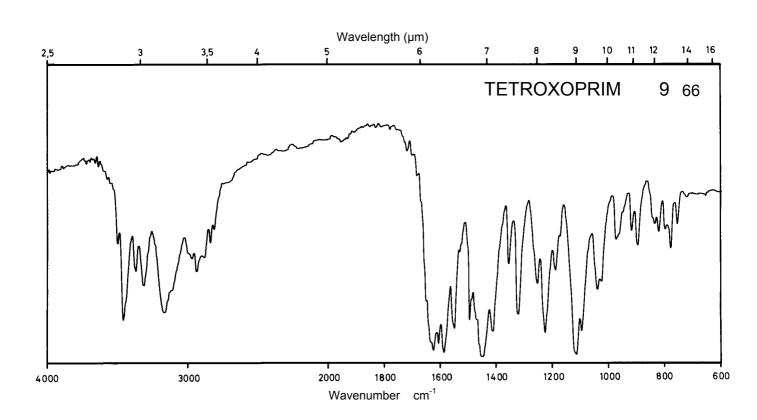





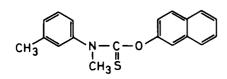


 M_r 307.4

Concentration 1 mg / 100 ml 10 mg / 100 ml

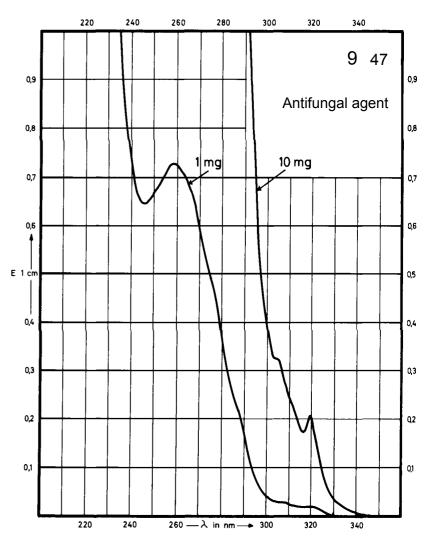

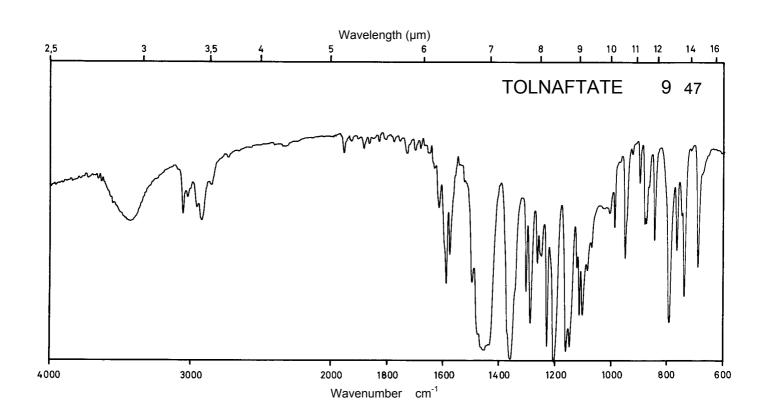





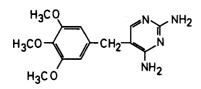


290.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm		271 nm	288 nm
E ^{1%} 1cm	232		208	252
ε	6730		6040	7320

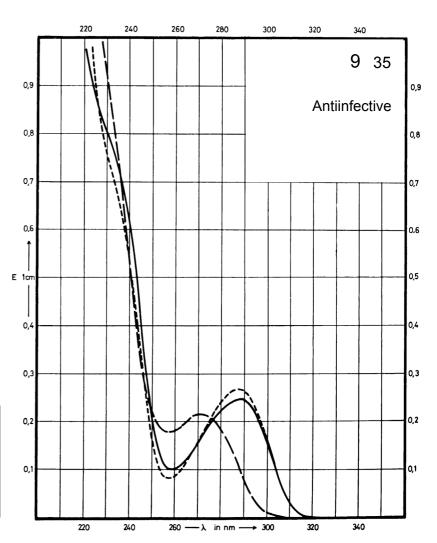

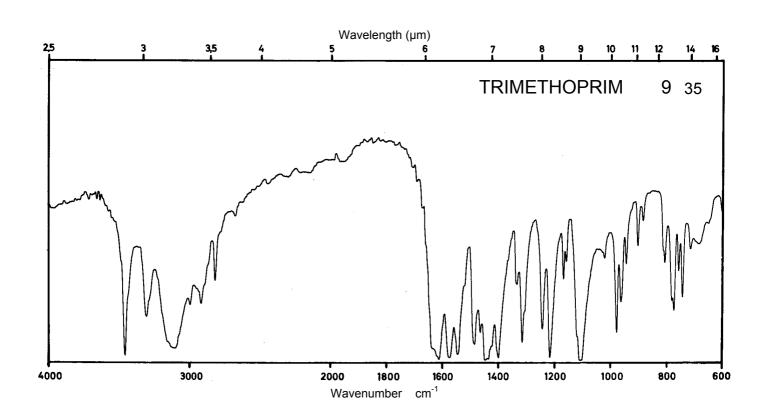





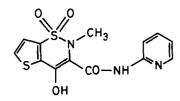


 M_r 337.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption			363 nm 263 nm	370 nm 257 nm
E ^{1%} 1cm			651 266	459 393
ε			22000 9000	15500 13250

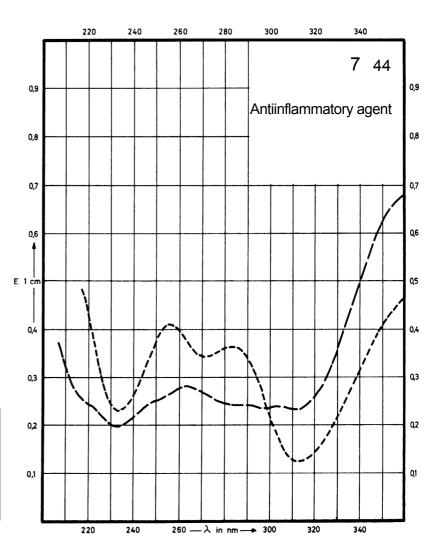

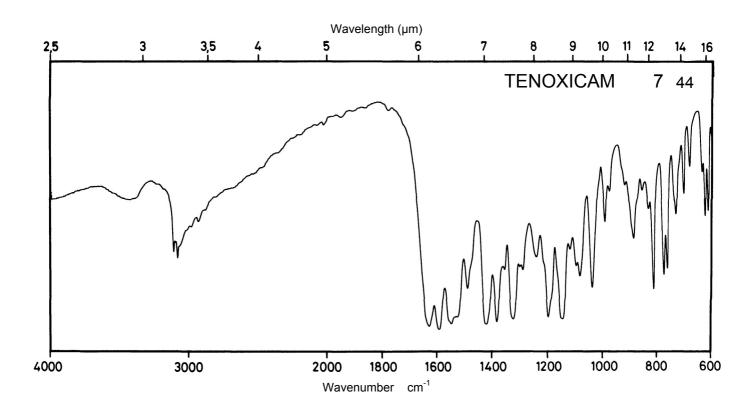



Name TIAPROFENIC ACID

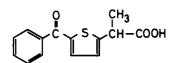


 M_r 260.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	303 nm 262 nm	315 nm 267 nm	308 nm	315 nm 266 nm
E ^{1%} 1cm	570 365	602 330	595	608 322
ε	14840 9500	15670 8590	15490	15830 8380

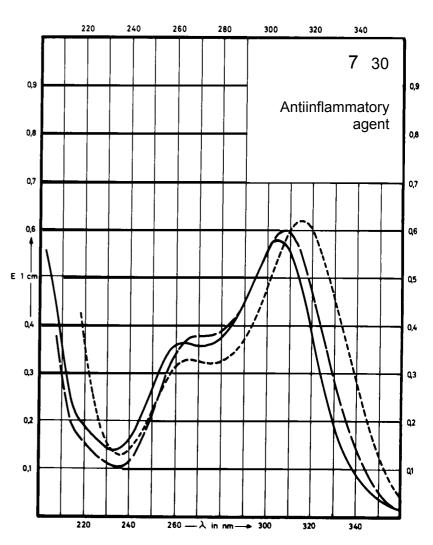

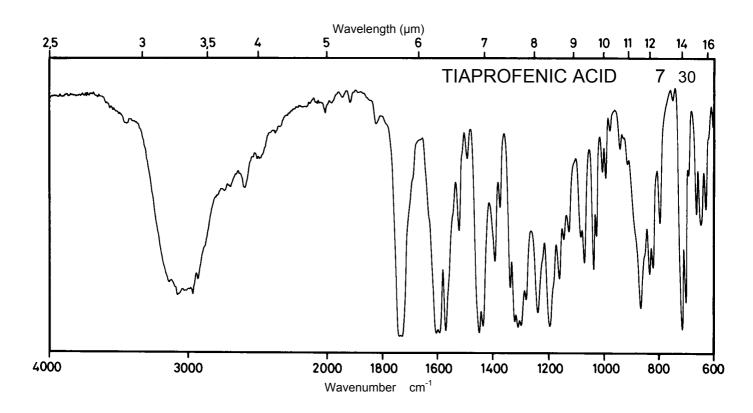





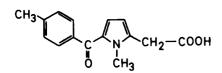


257.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm 255 nm		315 nm 262 nm	323 nm 260 nm
E ^{1%} 1cm	Decom- position observed		735 345	762 333
ε			18910 8880	19610 8570

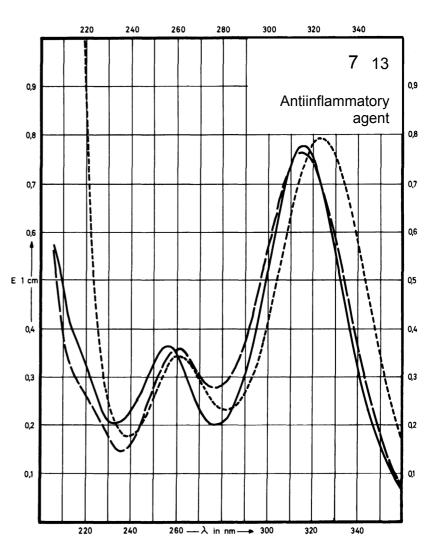

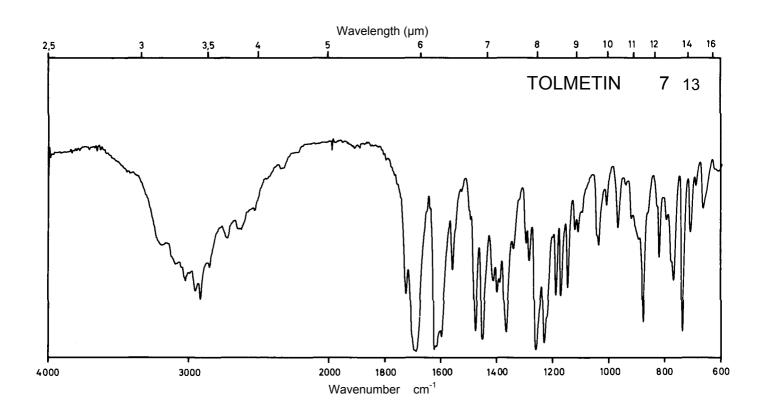




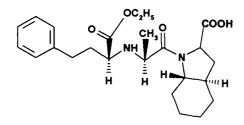


 M_r 430.6


Concentration 1 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm 258 nm 253 nm			
E ^{1%} 1cm	3.1 5.0 4.8			
ε	134 215 207			

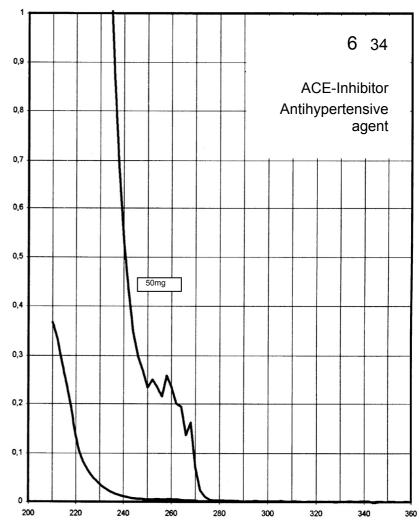

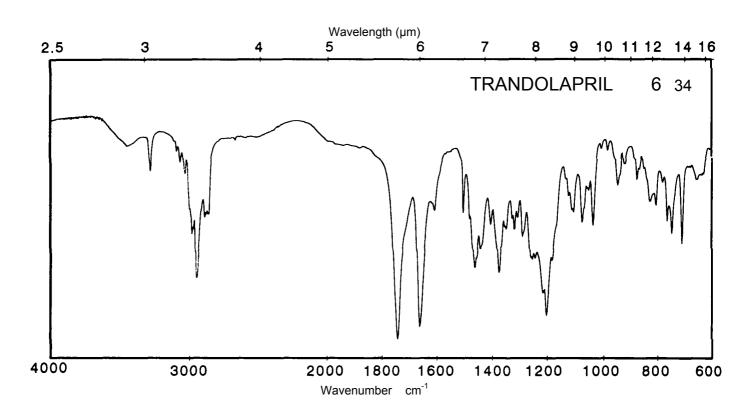



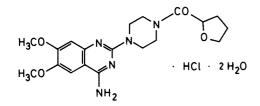



459.9

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	343 nm 247 nm		330 nm 245 nm	343 nm 249 nm
E ^{1%} 1cm	178 1230		233 1200	125 1310
ε	8200 56500		10700 55200	5740 60400

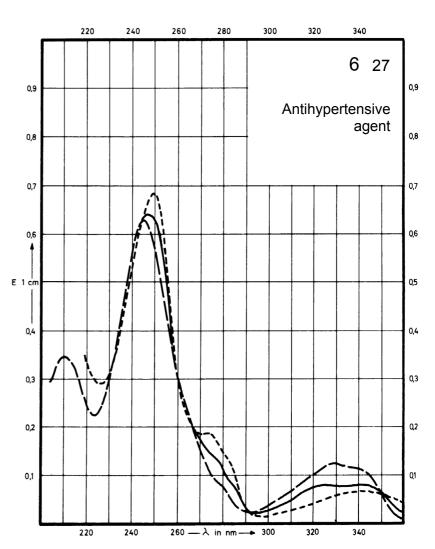

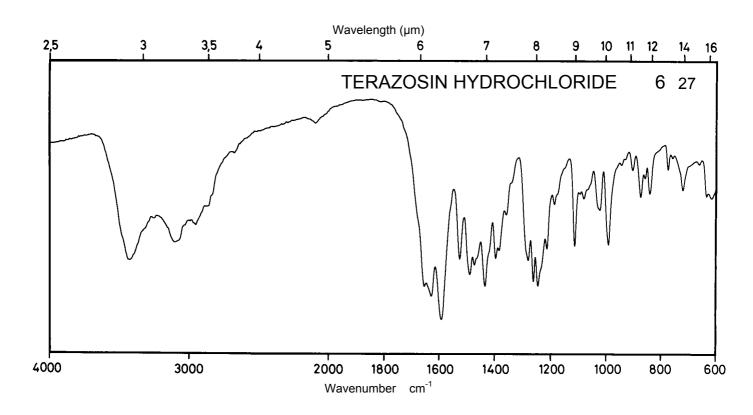


Name TIAMENIDINE HYDROCHLORIDE

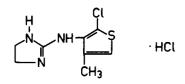



 M_r


252.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	242 nm	241 nm	241 nm	
E ^{1%} 1cm	297	296	296	
ε	7500	7500	7500	

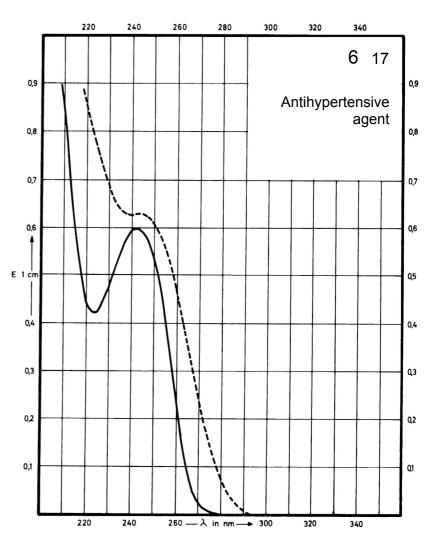

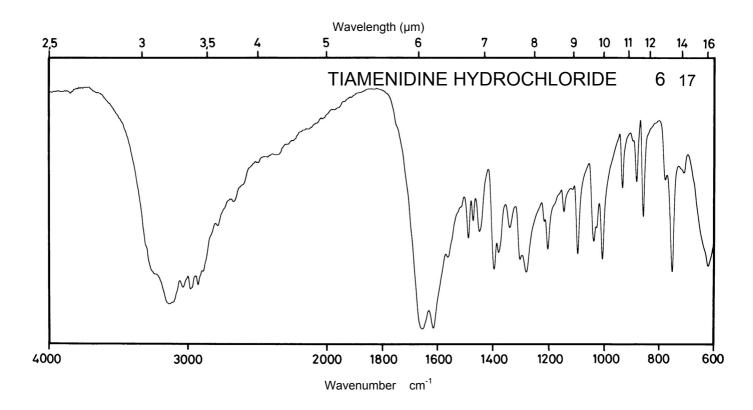

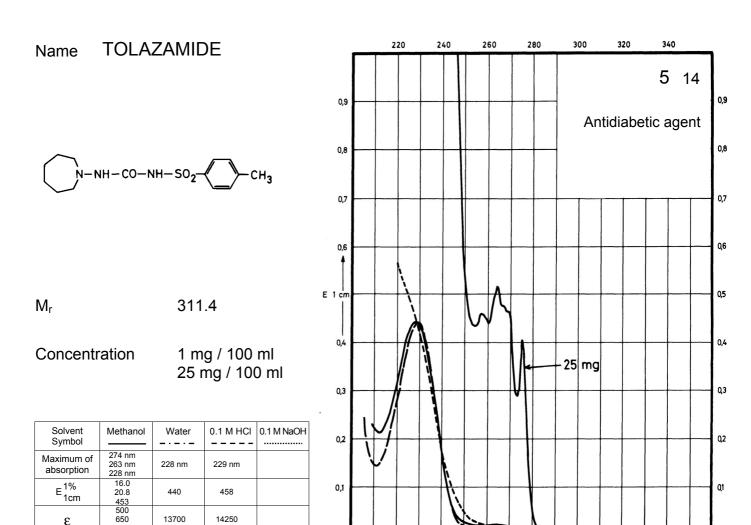

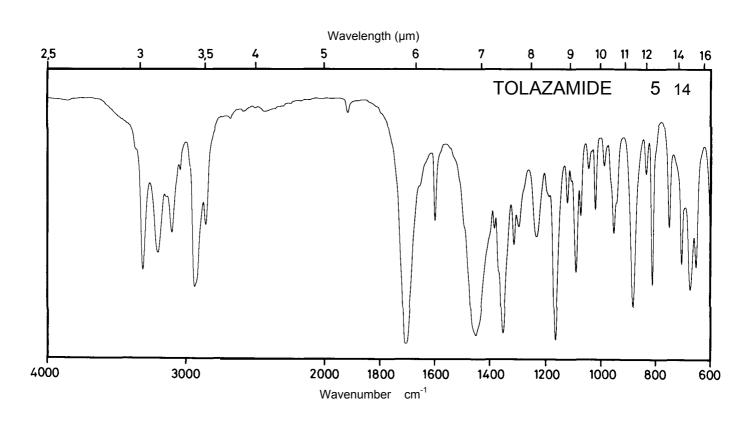

220

240

260 -

—λ in nm — 300


340

ε

14100

14250




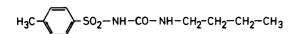


 M_r 270.4


Concentration 1 mg / 100 ml 25 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 228 nm		228 nm	
E ^{1%} 1cm	22.2 496		520	
ε	600 13410		14060	

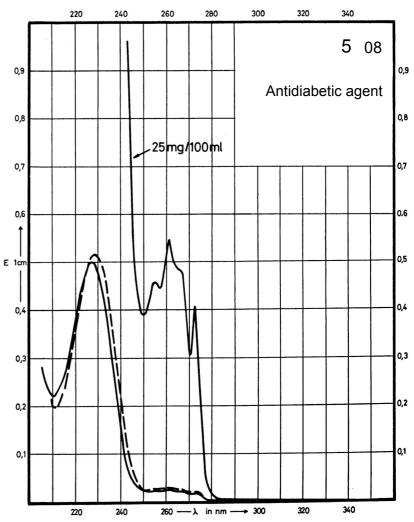

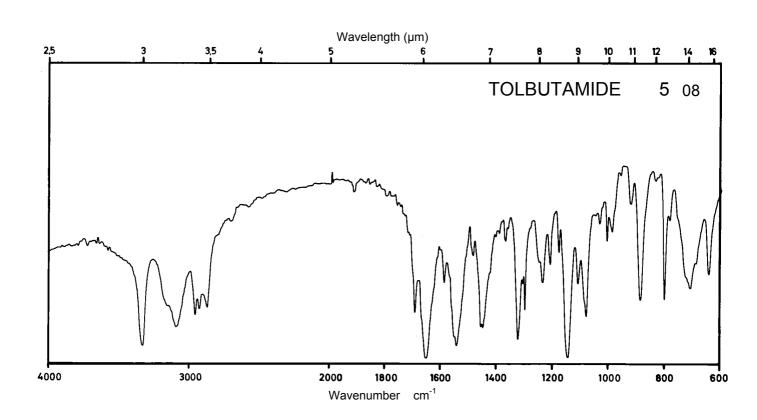



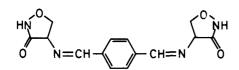



302.3

Concentration


0.76 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm		262 nm	270 nm
E ^{1%} 1cm	Decom- position observed		704	600
ε			21270	18150

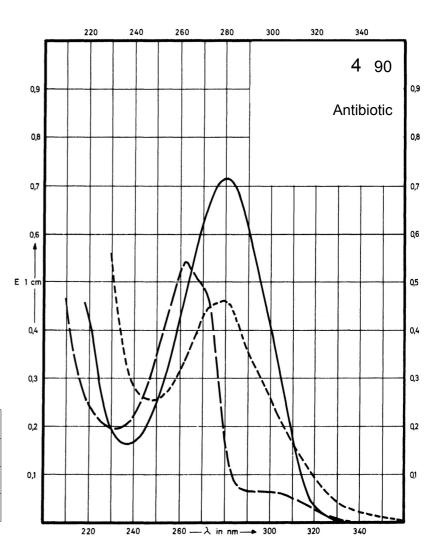

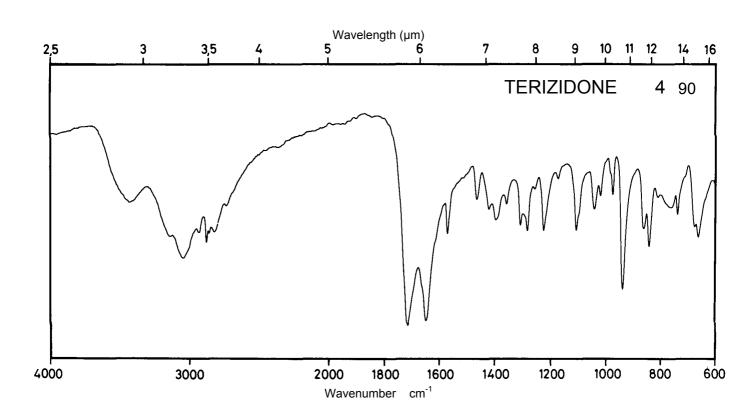




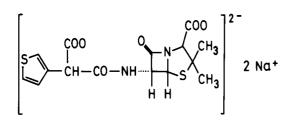


428.4


Concentration 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	232 nm		231 nm	233 nm
E ^{1%} 1cm	143		Decom- position observed	136
ε				

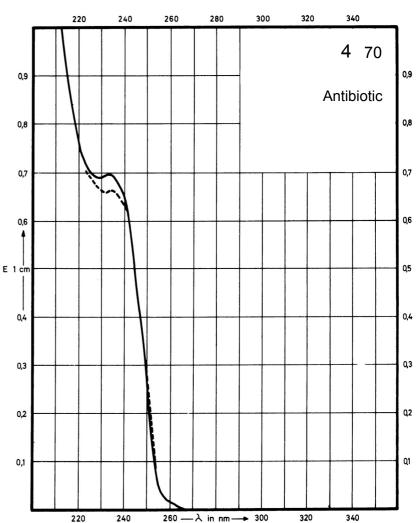

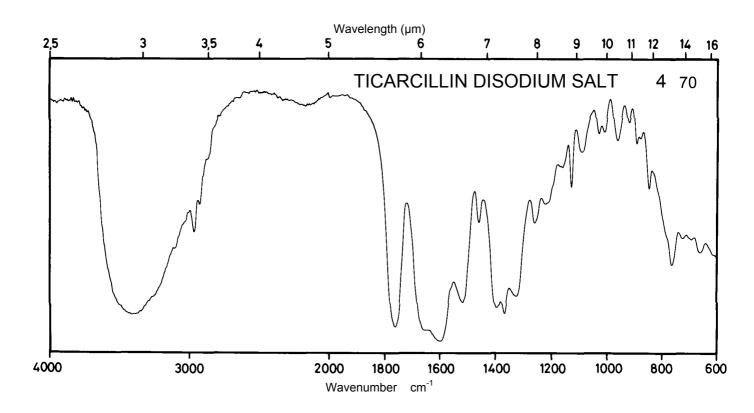

Name TETRACYCLINE

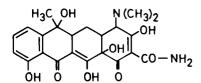



 M_r

Concentration 1.5 mg / 100 ml


444.4

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 363 nm 270 nm 356 nm 269 nm 380 nm 268 nm absorption E^{1%} 1cm 367 366 326 410 16300 16270 14490 18230 17670 14190 ε

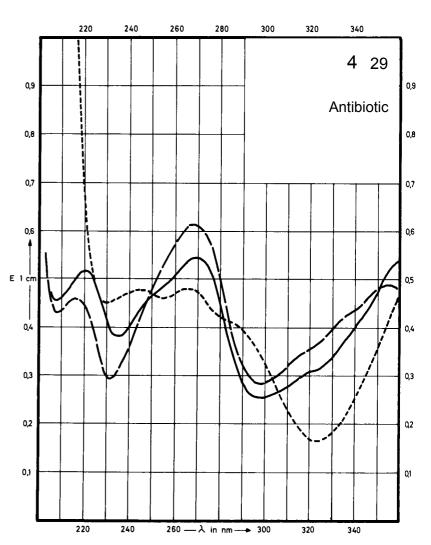

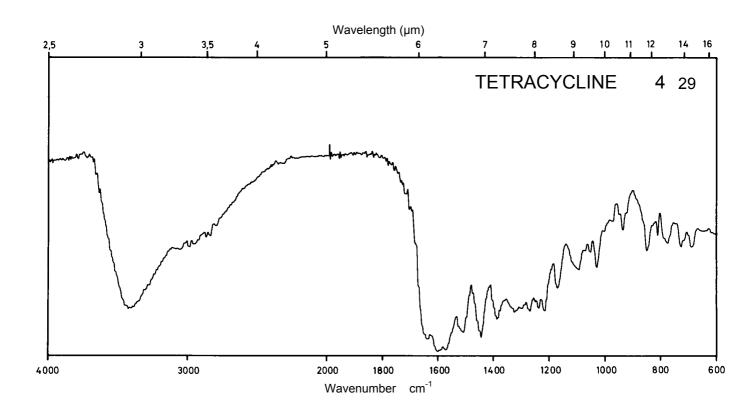

Name TETRACYCLINE HYDROCHLORIDE

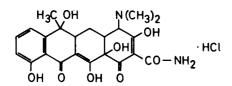


 M_{r}


480.9

Concentration

1.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	363 nm 270 nm		356 nm 269 nm	380 nm 268 nm
E ^{1%} 1cm	331 345		301 382	366 299
ε	15940 16590		14480 18380	17580 14360

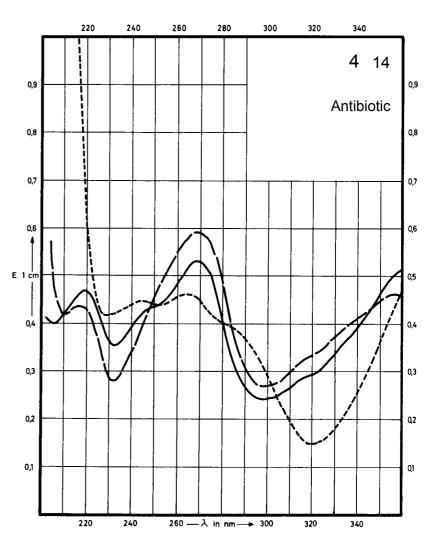

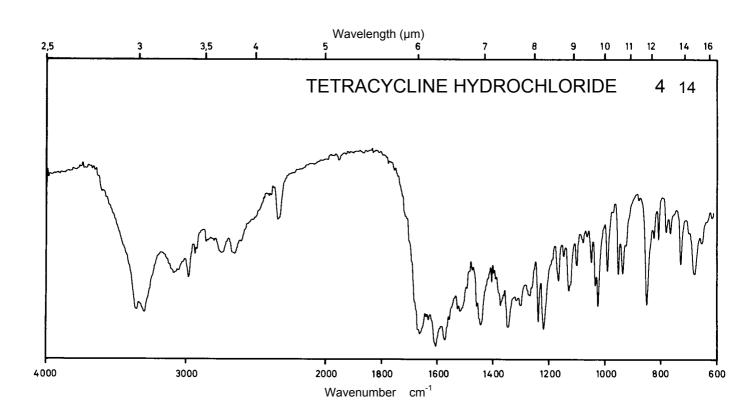


Name THIETHYLPERAZINE DIMALEATE


 M_{r}

631.8


Concentration 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 263 nm		311 nm 260 nm	
E ^{1%} 1cm	73 592		68 555	
3	4600 37400		4300 35100	

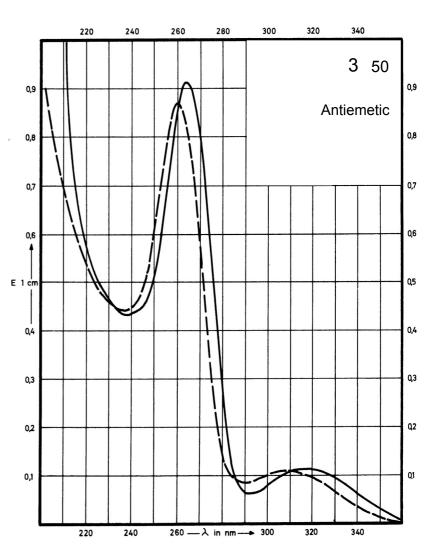

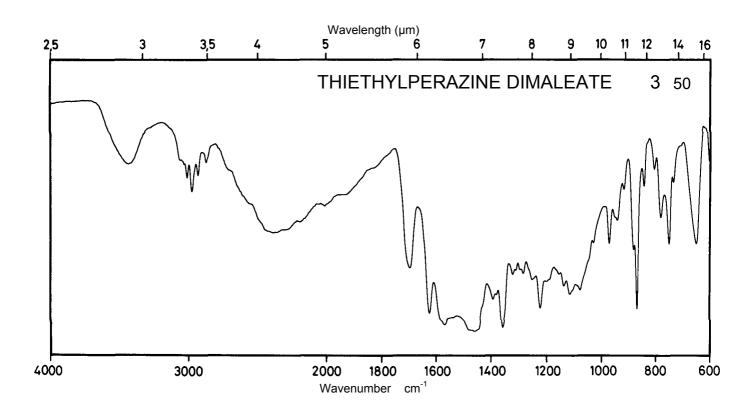



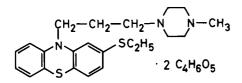



667.8

Concentration


1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	316 nm 263 nm		311 nm 260 nm	
E ^{1%} 1cm	68 547		64 518	
ε	4600 36500		4300 34600	

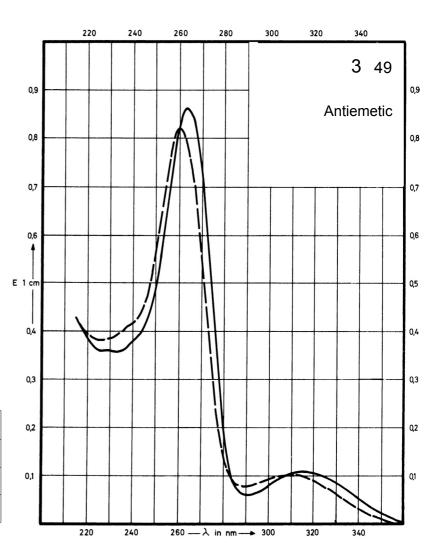

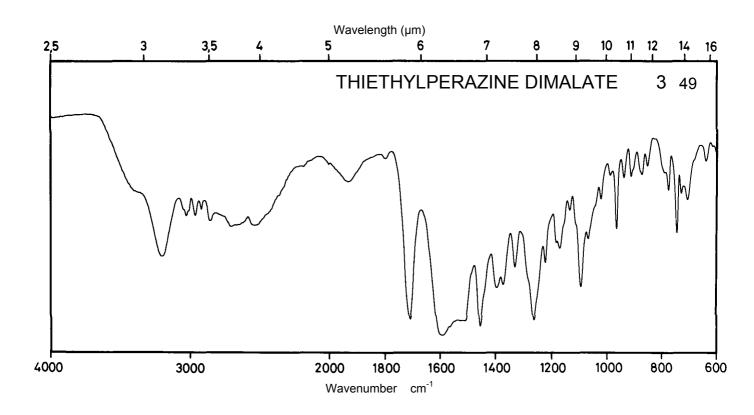



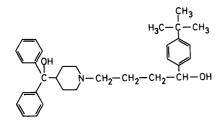



 M_r

471.7

Concentration 5 mg / 100 ml 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	259 nm			
E ^{1%} 1cm	14.3			
ε	670			

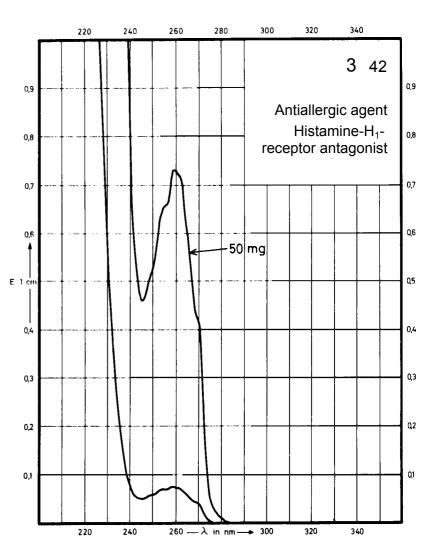

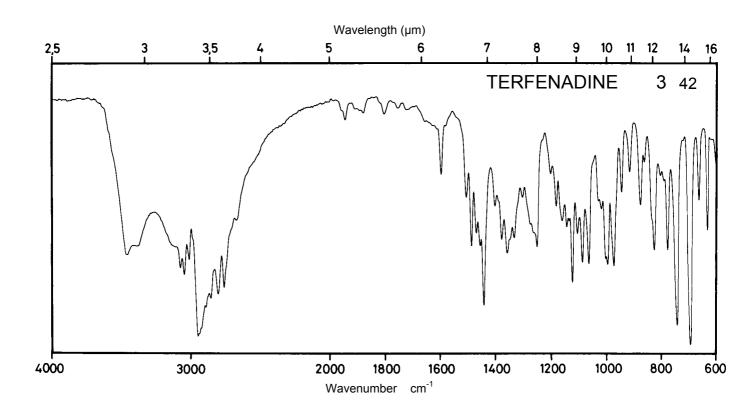

Name TRIPROLIDINE HYDROCHLORIDE

$$CH_3$$
 $C=CH-CH_2-N$ \cdot $HC1$ \cdot H_2O


 M_{r}


332.9

Concentration


1.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	280 nm 233 nm		290 nm	278 nm 228 nm
E ^{1%} 1cm	240 480		288	218 485
ε	7990 15980		9600	7260 16160

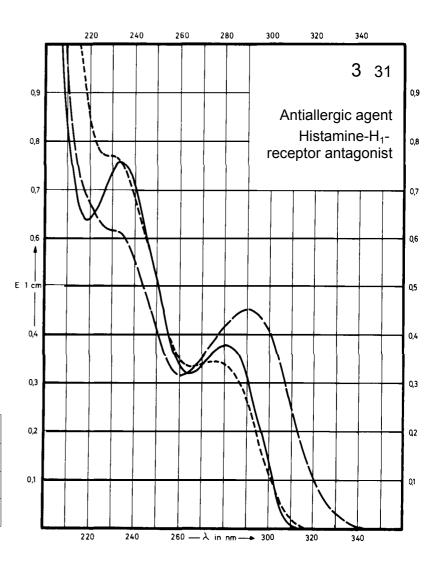

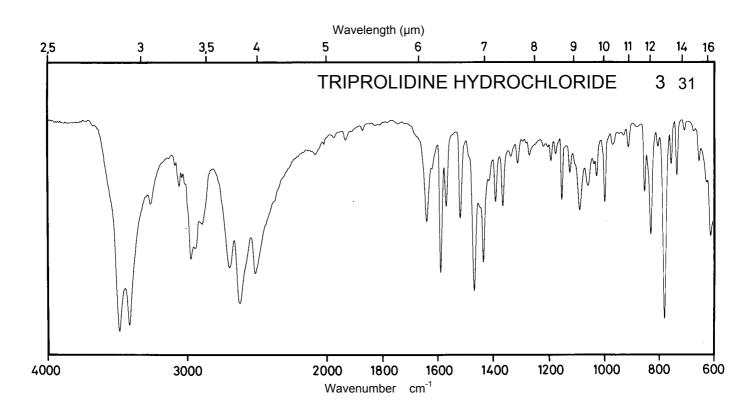



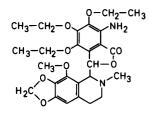



500.6

Concentration 4 m

4 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	340 nm			
E 1% 1cm	100			
ε	5000			

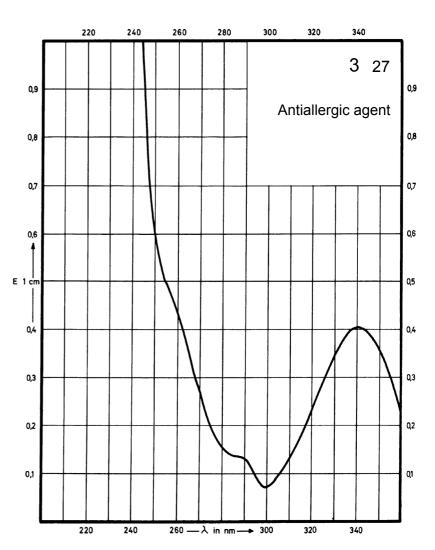

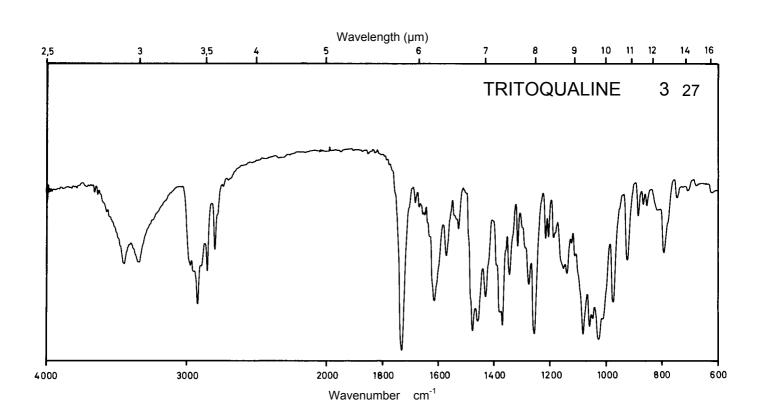


Name TRIPELENAMINE HYDROCHLORIDE


 M_r

291.8

Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	306 nm 245 nm		314 nm 239 nm	312 nm 249 nm
E ^{1%} 1cm	161 548		280 488	160 544
3	4700 16000		8150 14250	4670 15870

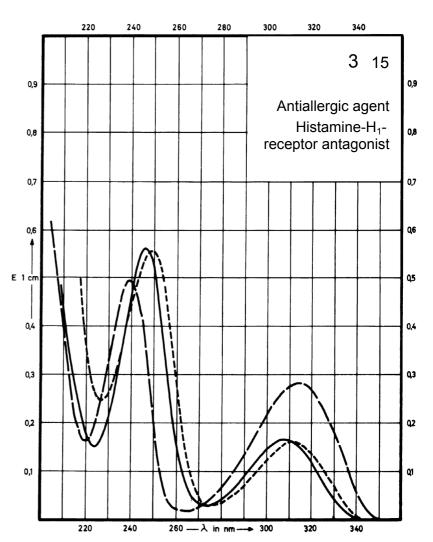

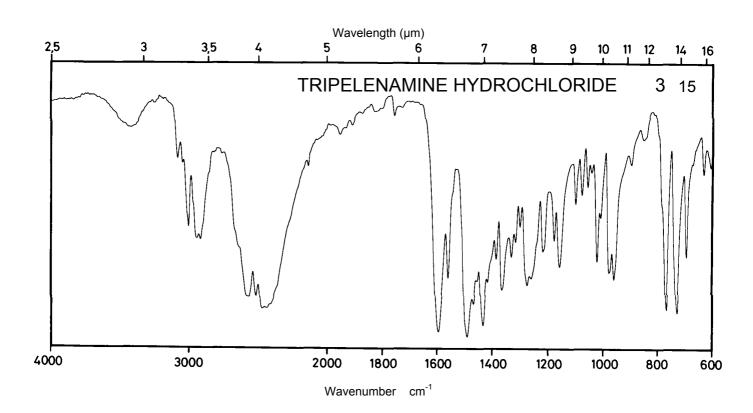


Name TRICLABENDAZOLE


 M_{r}

359.7


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	306 nm 254 nm		303 nm	
E ^{1%} 1cm	522 251		622	
3	18800 9000		22400	

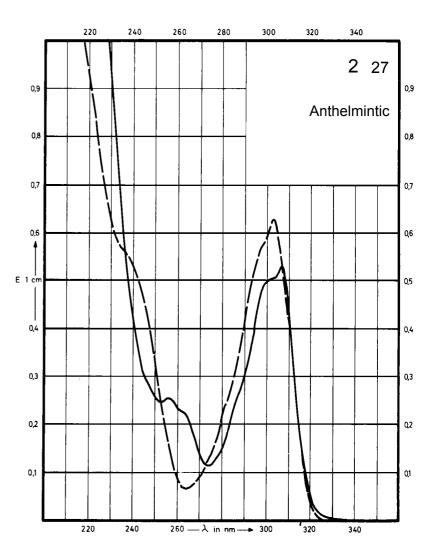

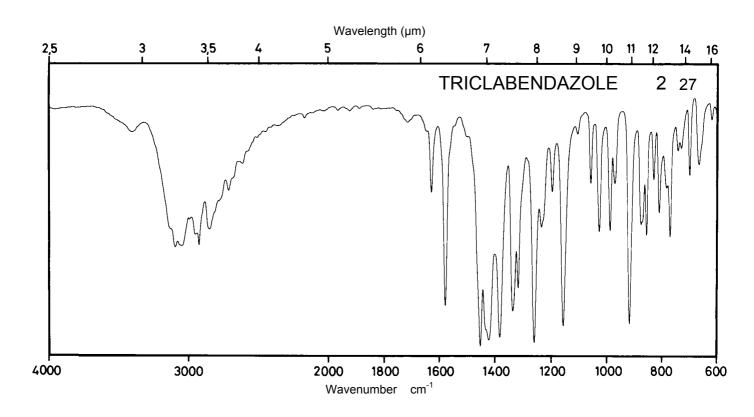



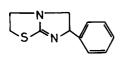



204.3

Concentration


1 mg / 100 ml 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm			
E ^{1%} 1cm	13			
3	265			

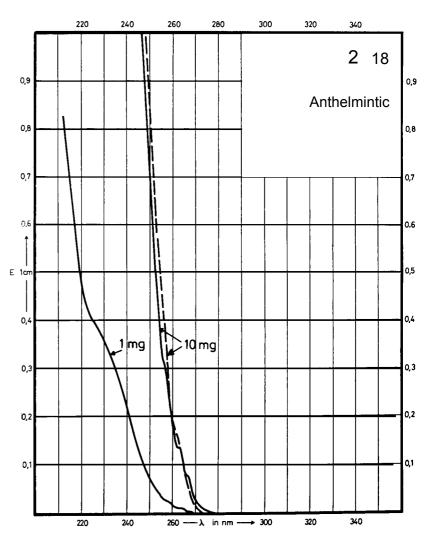

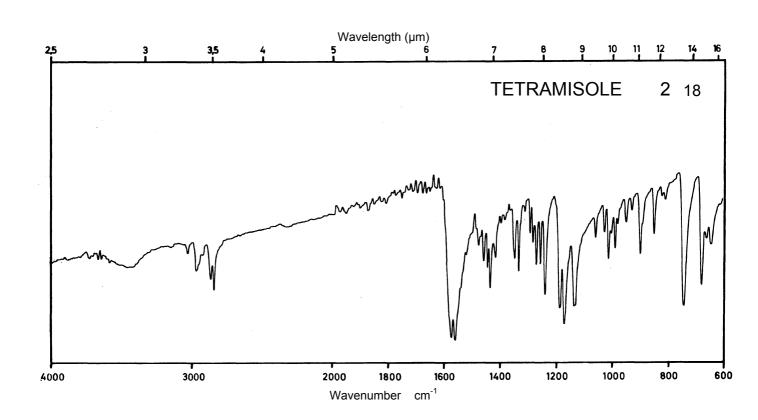

Name TIABENDAZOLE

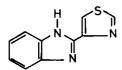


 M_{r}


201.3

Concentration


0.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	298 nm 241 nm		301 nm 242 nm	302 nm
E ^{1%} 1cm	1068 606		1253 637	1035
ε	21500 12200		25220 12820	20830

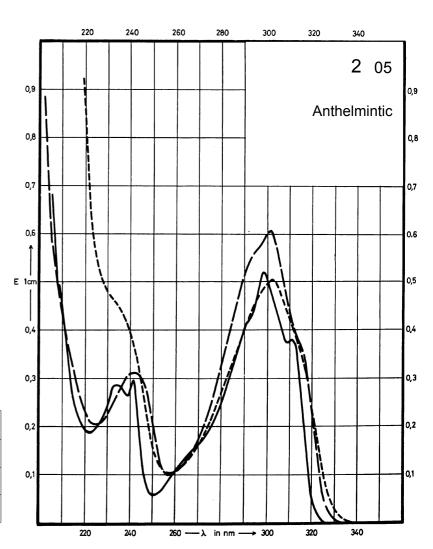

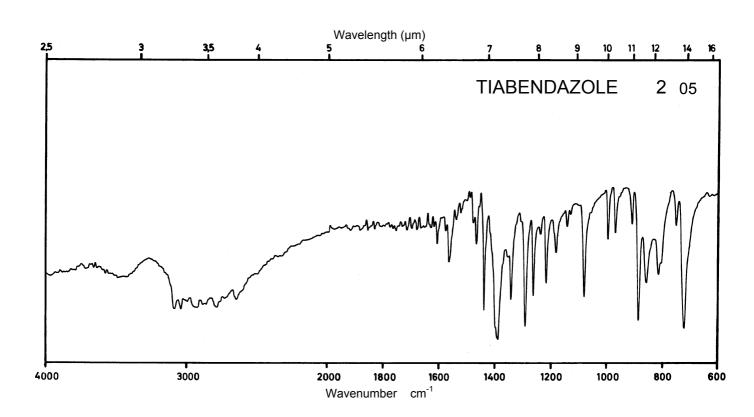


Name TRAMADOL HYDROCHLORIDE

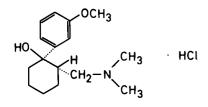


 M_r 299.8


Concentration 10 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	279 nm 273 nm		272 nm	272 nm
E ^{1%} 1cm	59 66		59	58
ε	1780 1960		1780	1720

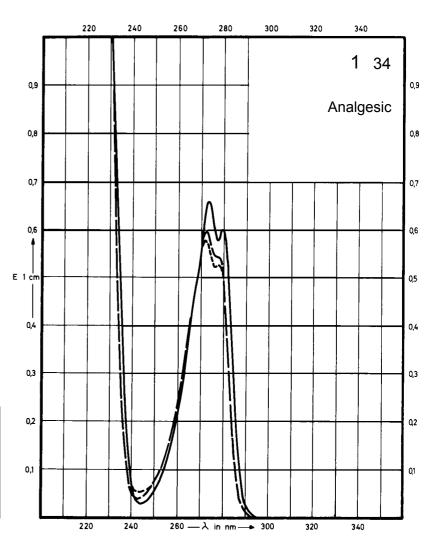

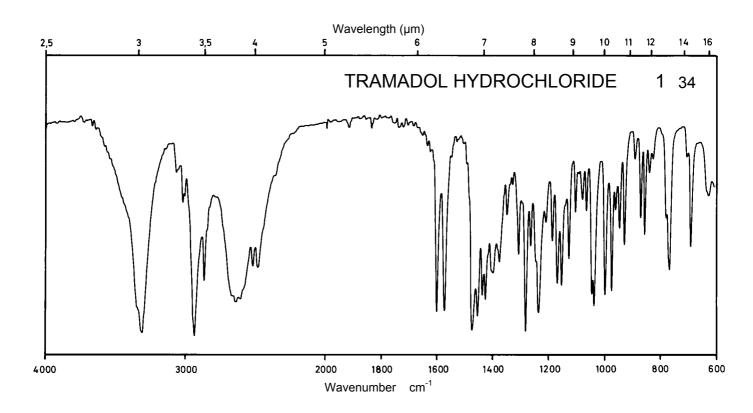



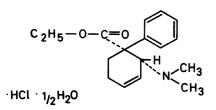

 M_r


318.8


Concentration 109 mg / 100 ml (CH₃OH)

101 mg / 100 ml


(HCI)

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm 257 nm 251 nm		262 nm 257 nm 250 nm	
E ^{1%} 1cm	5.0 6.3 5.2		5.2 6.5 5.3	
3	160 200 165		165 210	

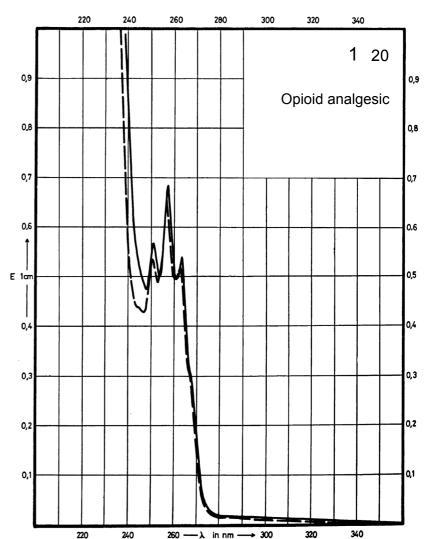

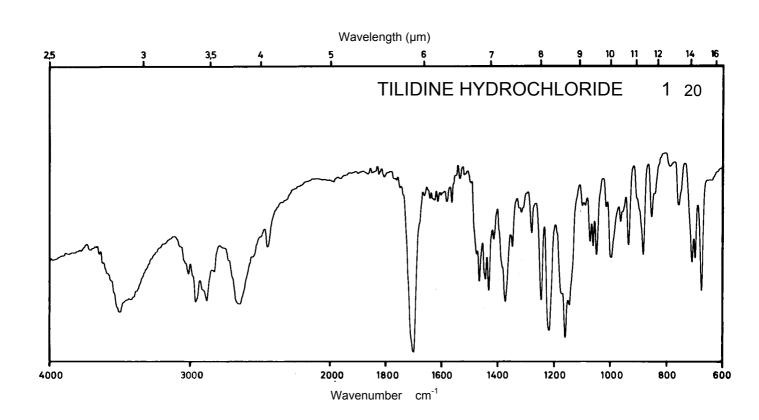




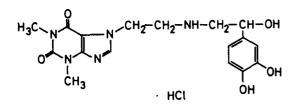



 M_r 411.8


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm		275 nm	Decom- position observed
E ^{1%} 1cm	272		277	
ε	11200		11400	

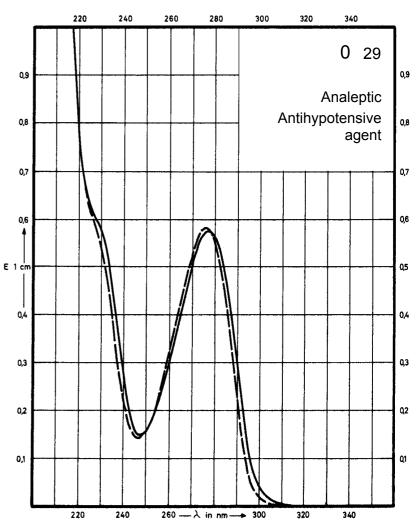


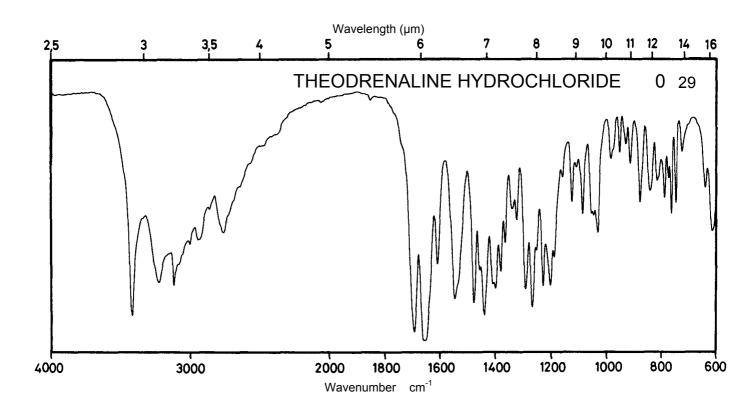

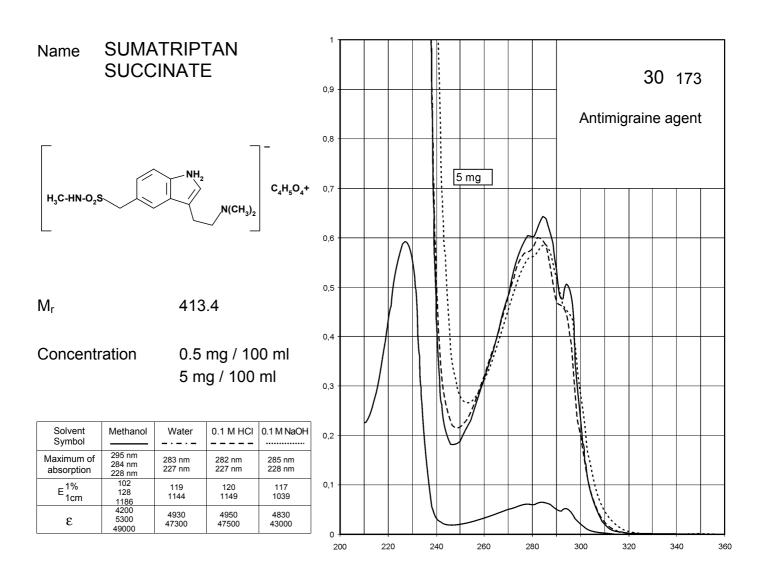

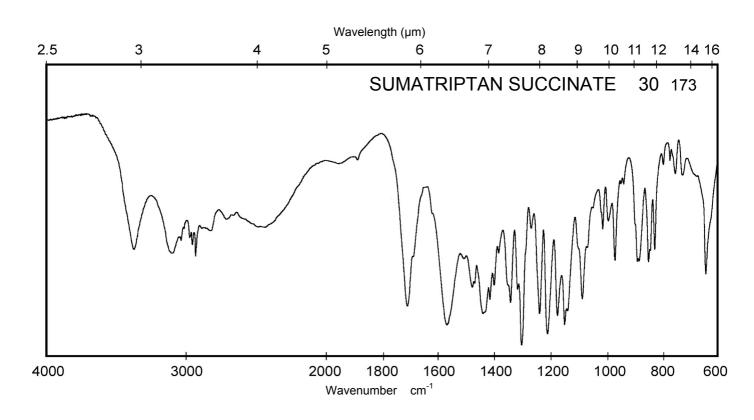




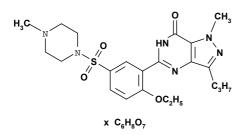



 M_{r} 666.8


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	291 nm	293 nm 225 nm	292 nm 227 nm	290 nm
E ^{1%} 1cm	217	199 455	197 462	169
ε	14500	13300 30300	13100 30800	11300

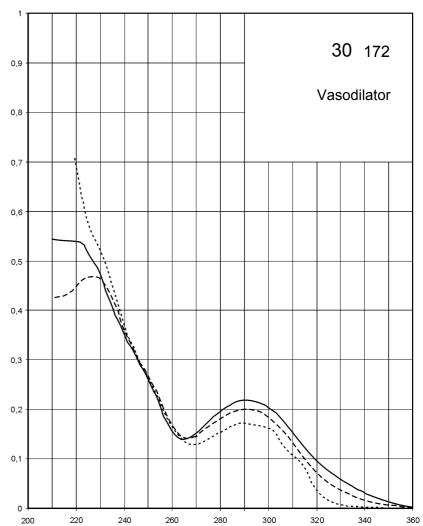

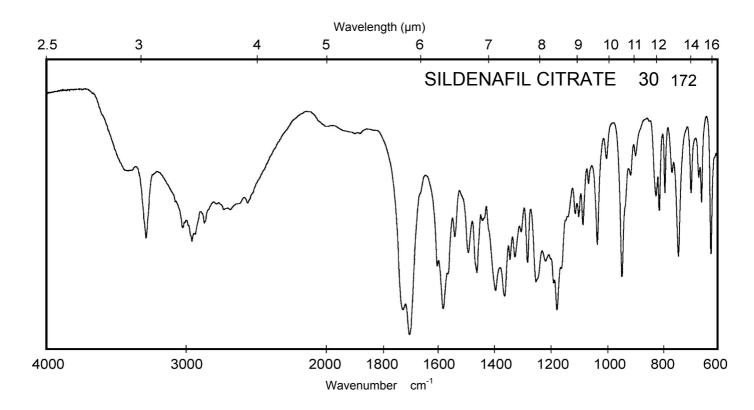

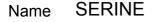


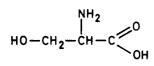



 M_r 105.1


Concentration 100 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

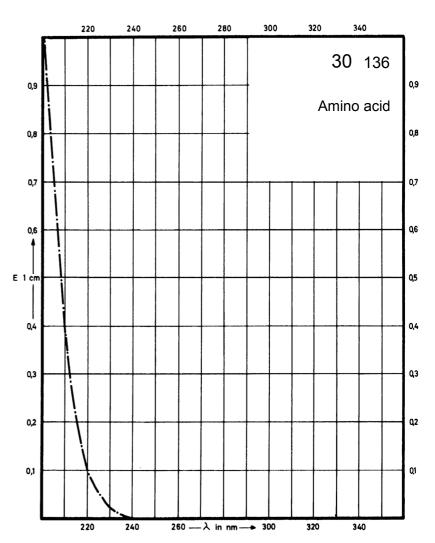

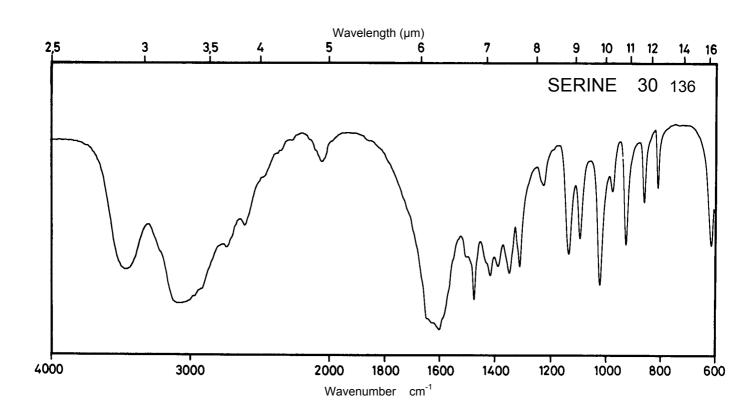

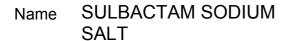


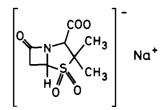


255.2

Concentration 12 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				Decom- position observed
E ^{1%} 1cm				
3				

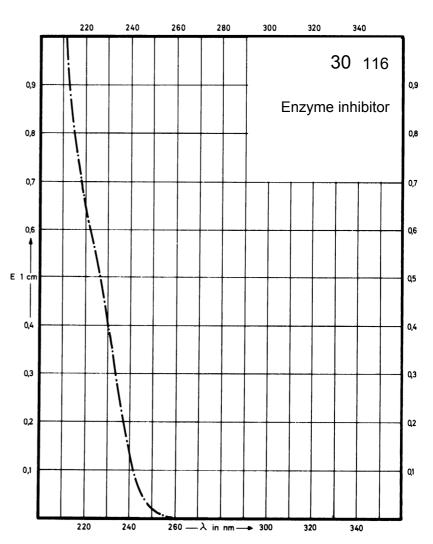

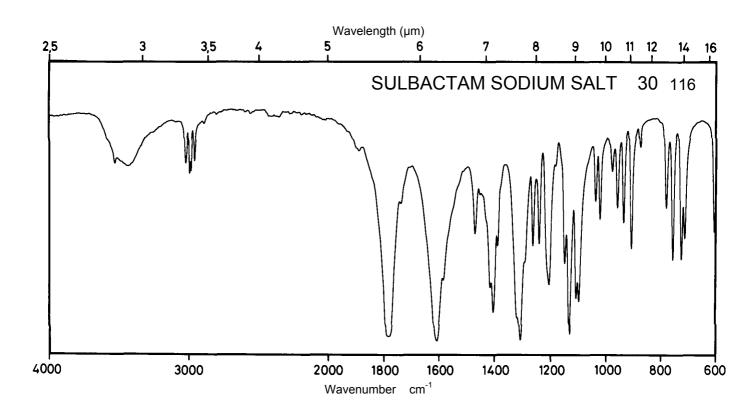

Name SOMATOSTATIN


Ala - Gly - Cys - Lys - Asn - Phe - Phe - Trp Cys - Ser - Thr - Phe - Thr - Lys


 M_{r}

1637.9


Concentration 1.8 mg / 100 ml 22 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		288 nm 280 nm		
E ^{1%} 1cm		25 30		
ε		4120 4870		

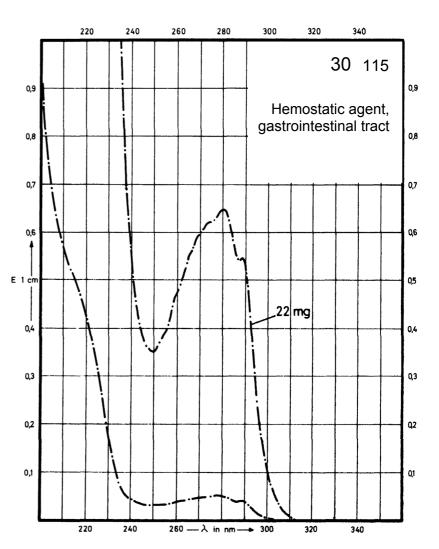

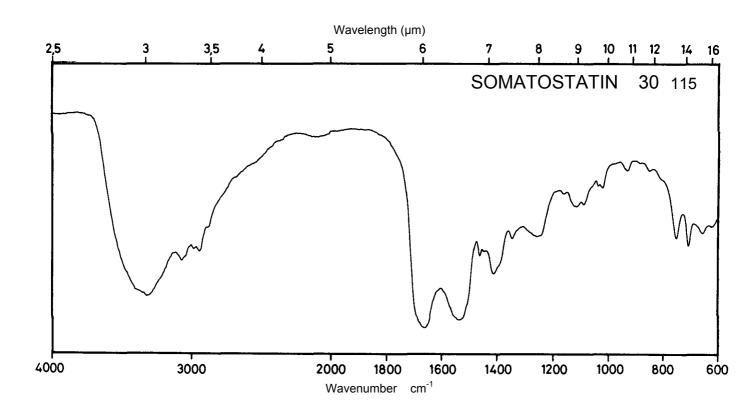




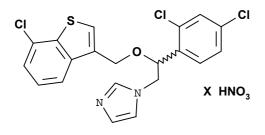


 M_r 500.8


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	302 nm 260 nm 225 nm		302 nm 260 nm 224 nm	
E ^{1%} 1cm	77 112 818		71 112 755	
ε	3850 5630 40900		3560 5610 37800	

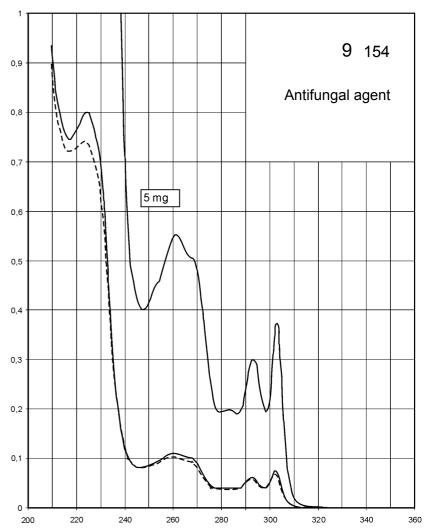

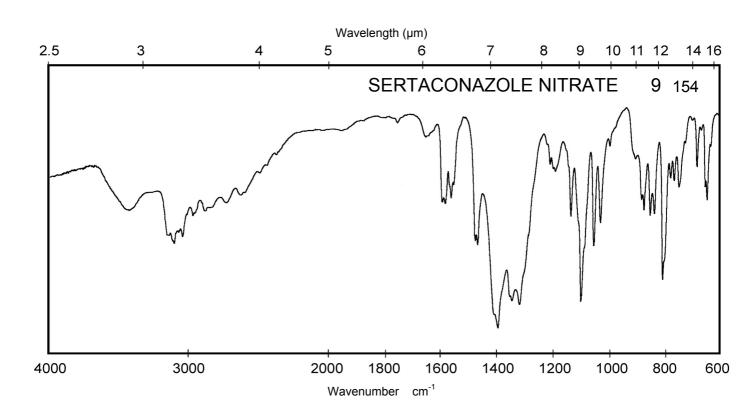




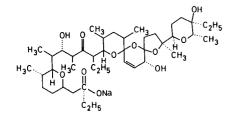


773.0


Concentration 50 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		287 nm		
E ^{1%} 1cm		0.65		
3		50		

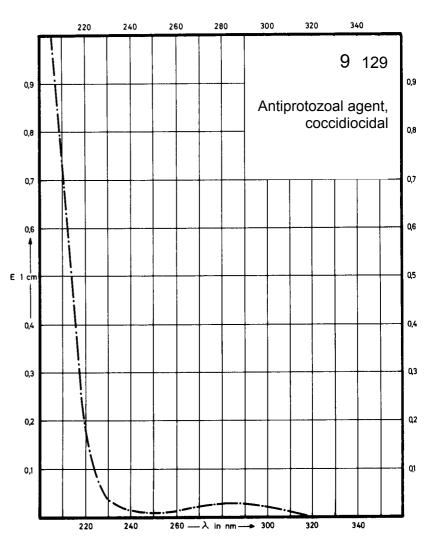

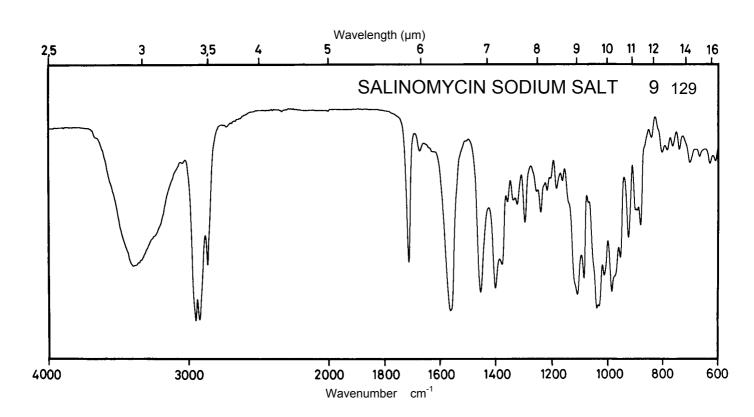




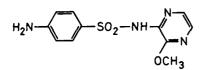


280.3


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm		304 nm 219 nm	315 nm 248 nm
E ^{1%} 1cm	730		292 535	324 795
3	20450		8180 15000	9100 22300

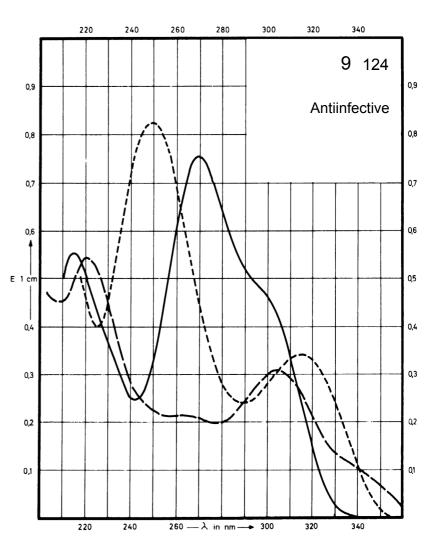

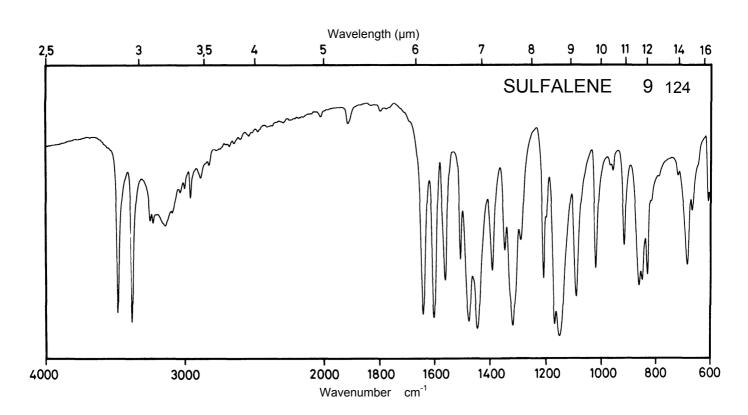




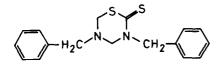


314.5


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	290 nm 248 nm			
E ^{1%} 1cm	381 275			
ε	12000 8650			

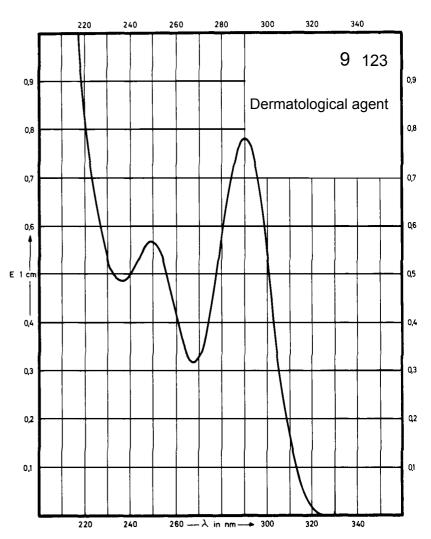

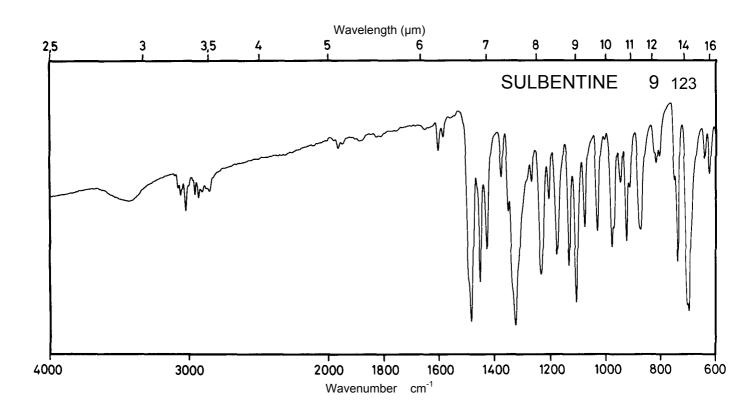



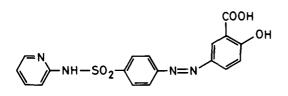



398.4

Concentration

0.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	368 nm 247 nm		350 nm 233 nm	456 nm 286 nm
E ^{1%} 1cm	680 393		662 540	750 359
ε	27110 15660		26400 21500	29860 14290

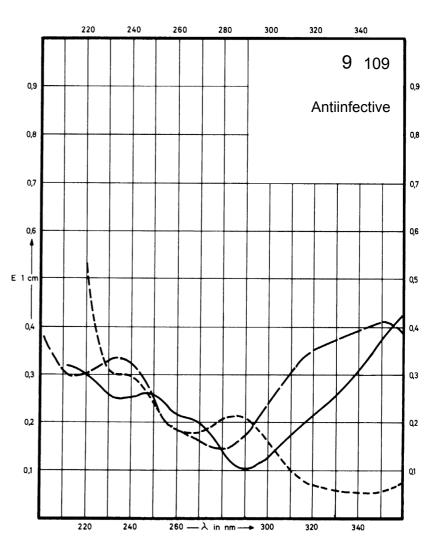

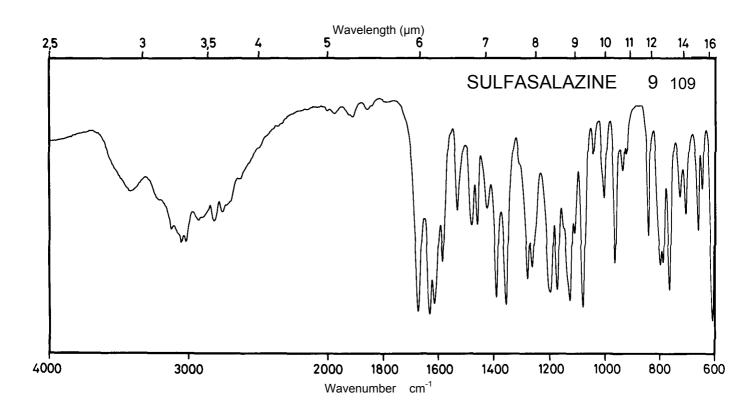





 $CH_3 - (CH_2)_{16} - COOH$


 M_r 284.5

Concentration 100 mg / 100 ml


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

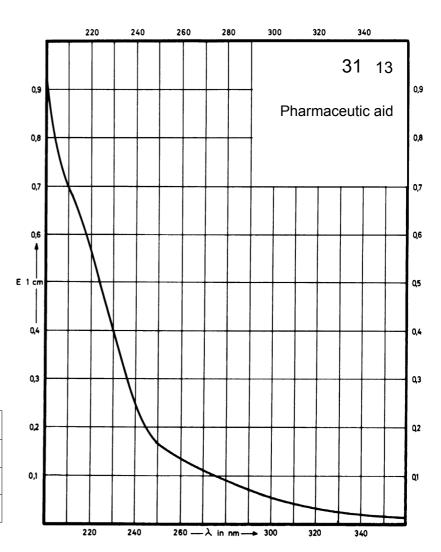


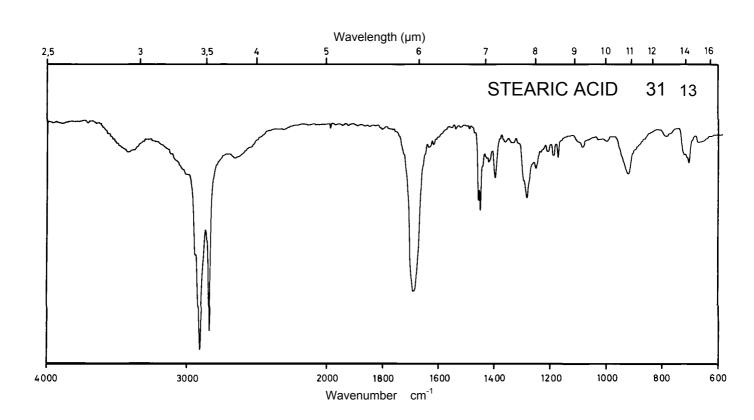

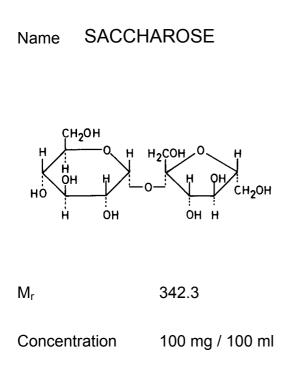


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

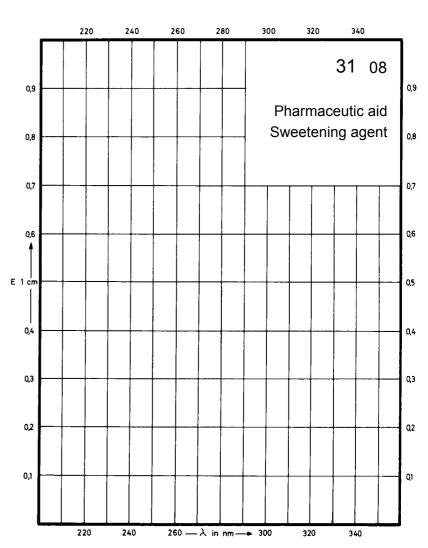

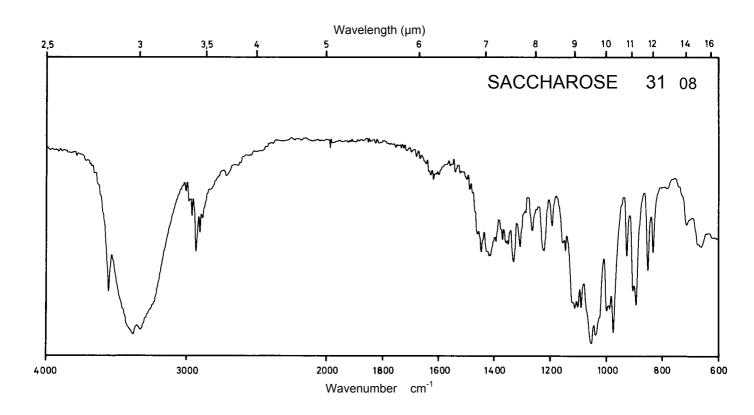




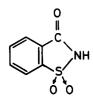


183.2

Concentration 8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm		278 nm	268 nm
E ^{1%} 1cm	71		60	73
ε	1300		1100	1340

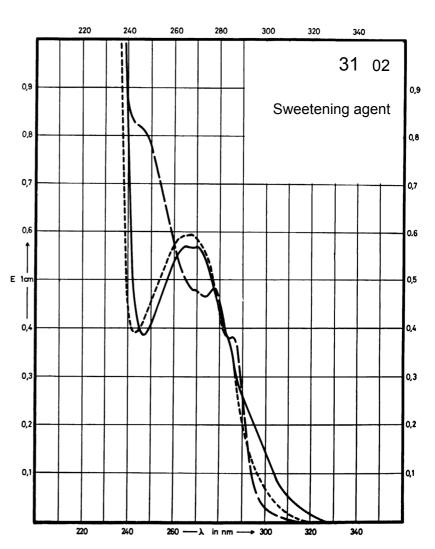

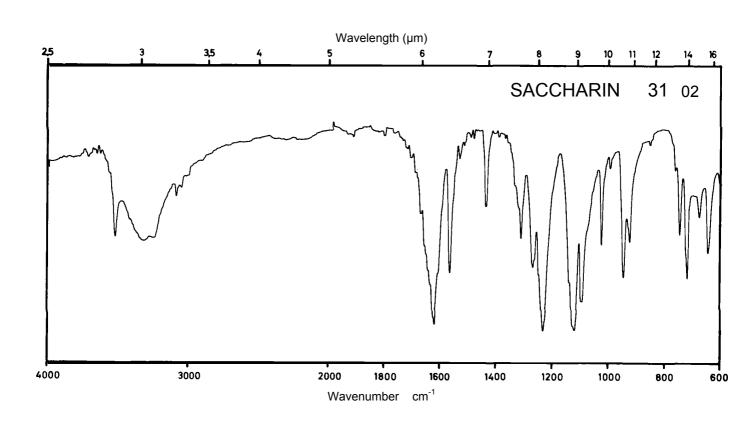



Name SULOCTIDIL

 M_r 337.6


Concentration 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm 219 nm		258 nm	
E ^{1%} 1cm	259 290		258	
ε	8750 9800		8700	

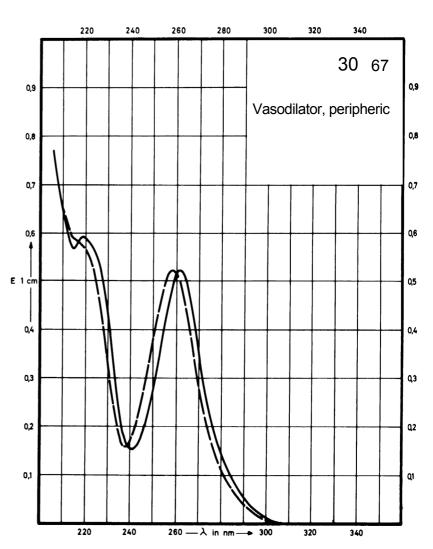

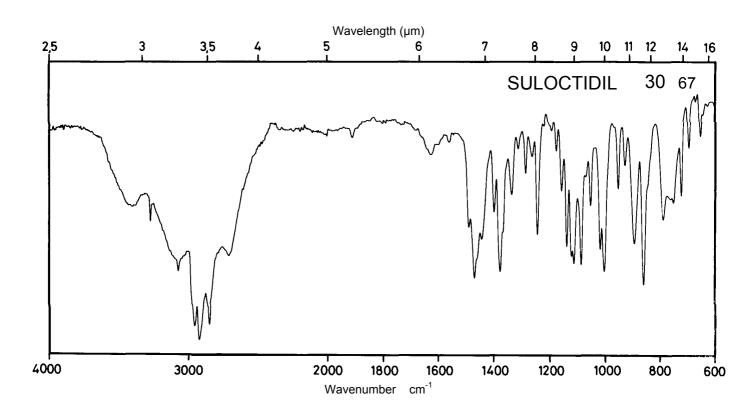



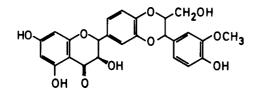



 M_r


Concentration 1 mg / 100 ml


482.4

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of absorption 287 nm 286 nm 310 nm E^{1%} 1cm 450 411 423 21700 19800 20400 ε

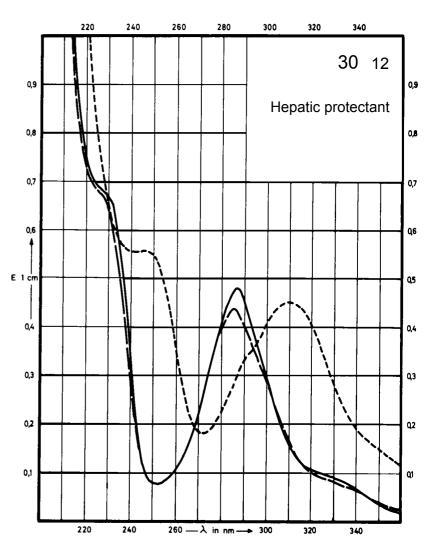


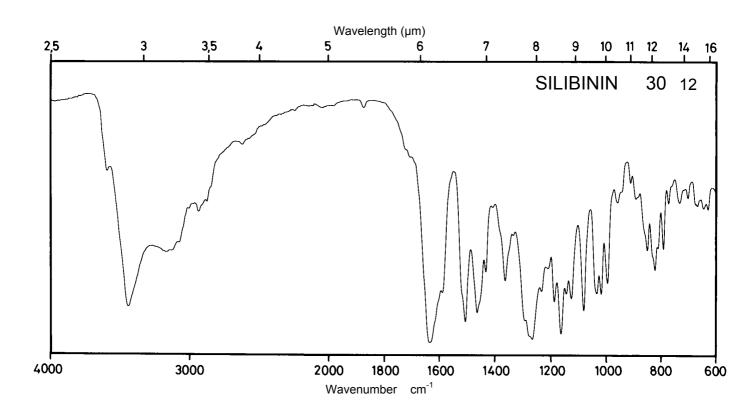

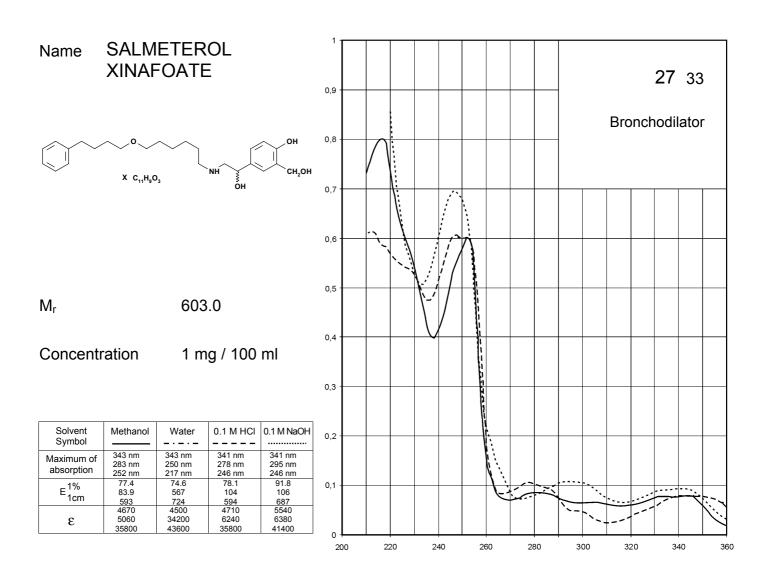

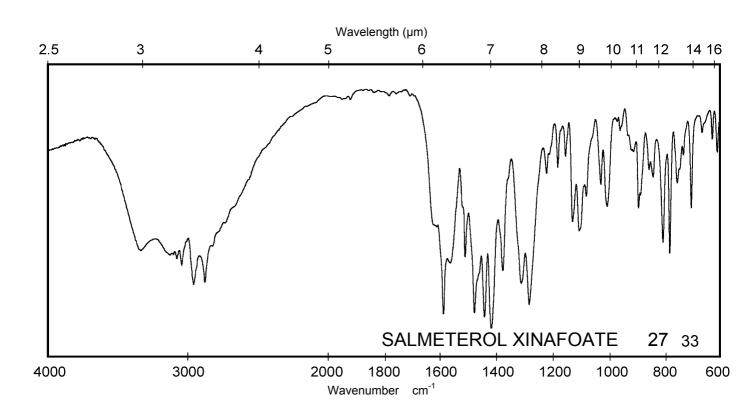




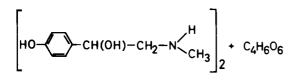


 M_r 484.5


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of 275 nm 224 nm 272 nm 222 nm 291 nm 241 nm absorption 67 412 100 559 3250 19960 2710 18170 4850 27080 ε

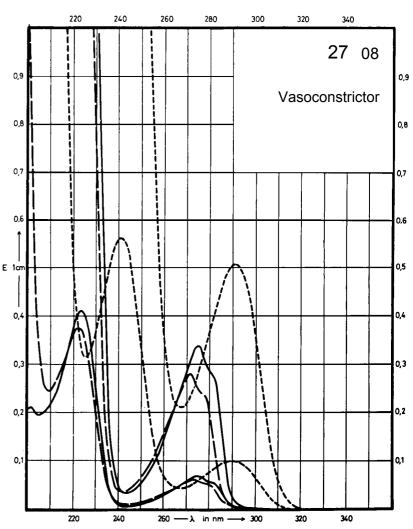

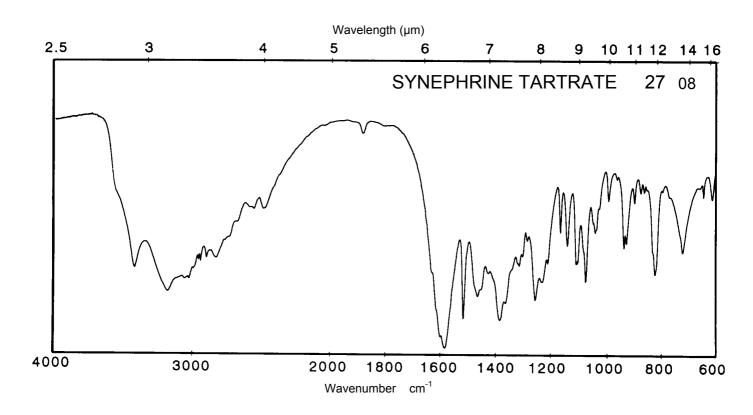




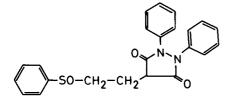


404.5


Concentration 1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm		237 nm	259 nm
E ^{1%} 1cm	502		447	574
ε	20330		18090	23210

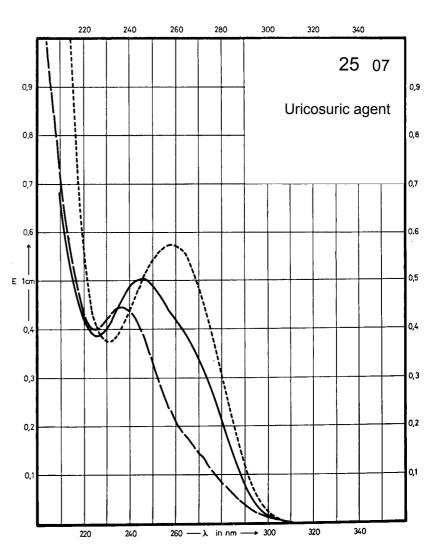

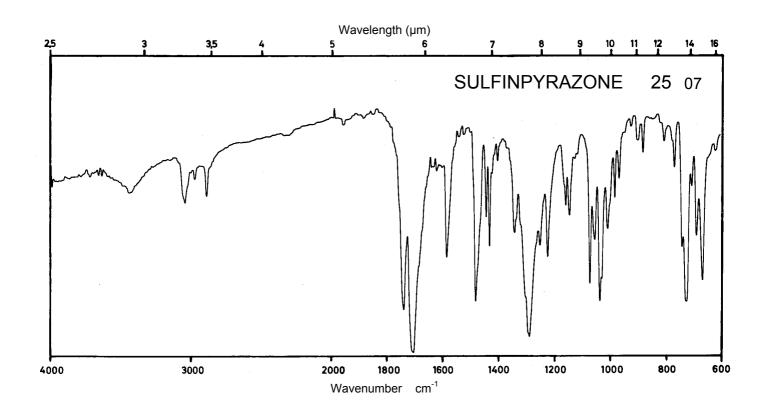


Name SPIRONOLACTONE

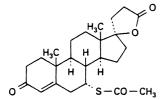


 M_{r}

416.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm		242 nm	292 nm 247 nm
E ^{1%} _{1cm}	471		464	118 428
3	19620		19330	4920 17830

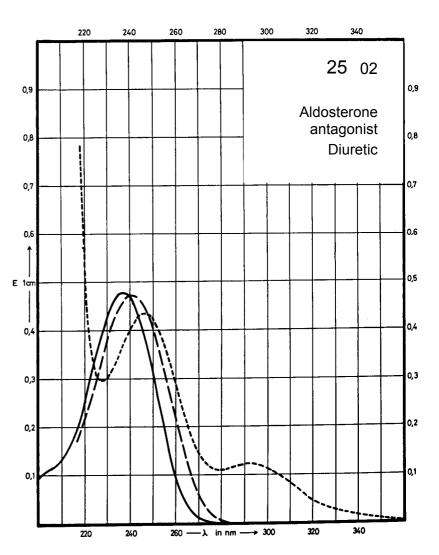

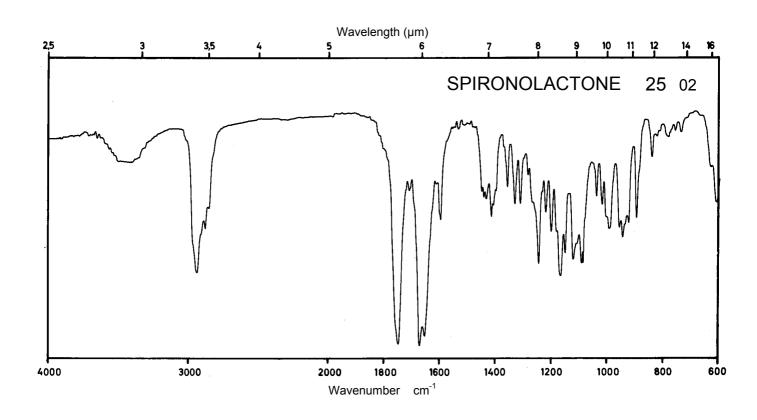



Name SALBUTAMOL SULFATE

 M_{r}

576.7

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	278 nm	276 nm	276 nm	295 nm
	227 nm	225 nm	225 nm	245 nm
E ^{1%}	68	61	61	110
1cm	289	272	272	423
ε	3930	3500	3500	6330
	16700	15700	15700	24400

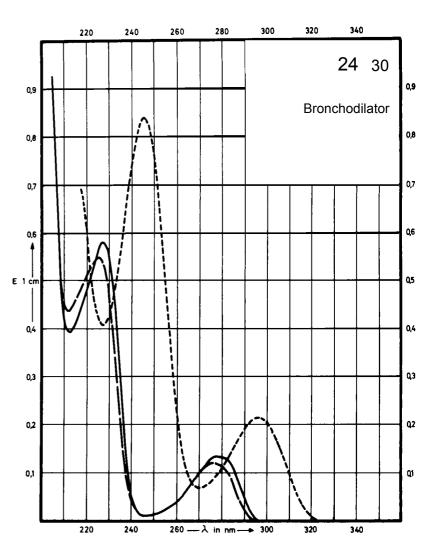

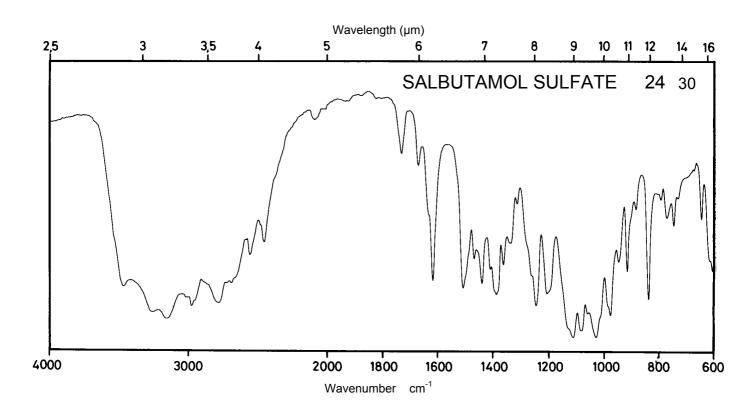



$$\begin{bmatrix} \mathsf{CH_2} - \mathsf{CO} - \mathsf{O} - \mathsf{CH_2} - \mathsf{CH_2} - \mathsf{N} \, (\mathsf{CH_3})_3 \\ \mathsf{CH_2} - \mathsf{CO} - \mathsf{O} - \mathsf{CH_2} - \mathsf{CH_2} - \mathsf{N} \, (\mathsf{CH_3})_3 \end{bmatrix}^{2+} \\ \quad \cdot \quad 2 \, \mathsf{H_2O}$$


 M_r


•

Concentration 100 mg / 100 ml


397.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

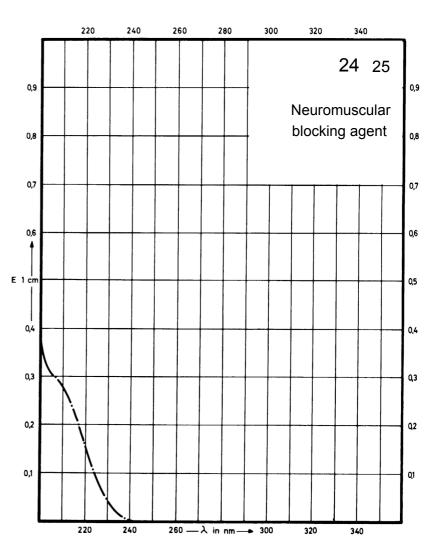

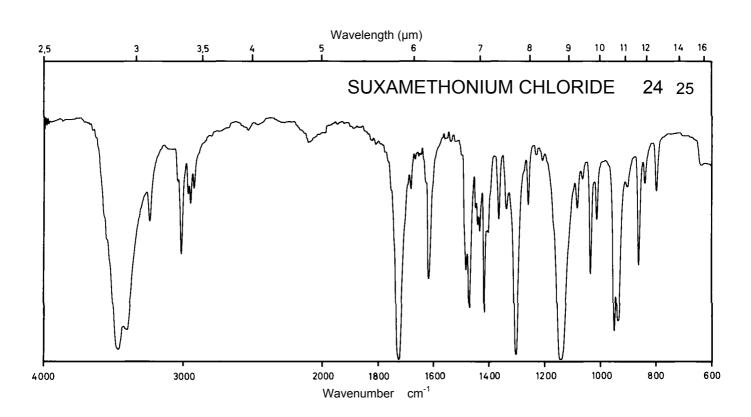





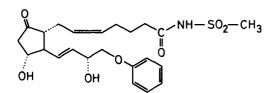


465.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	277 nm 270 nm			
E ^{1%} 1cm	62 73			
3	2900 3400			

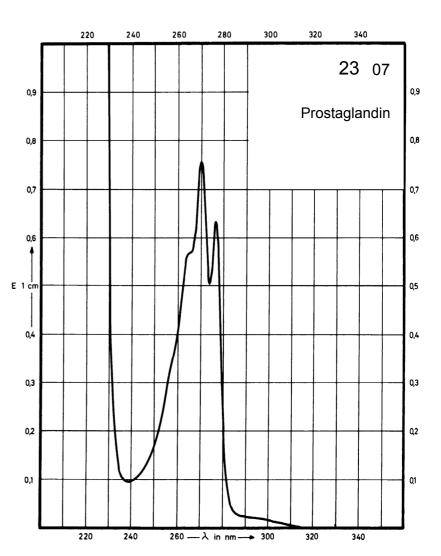

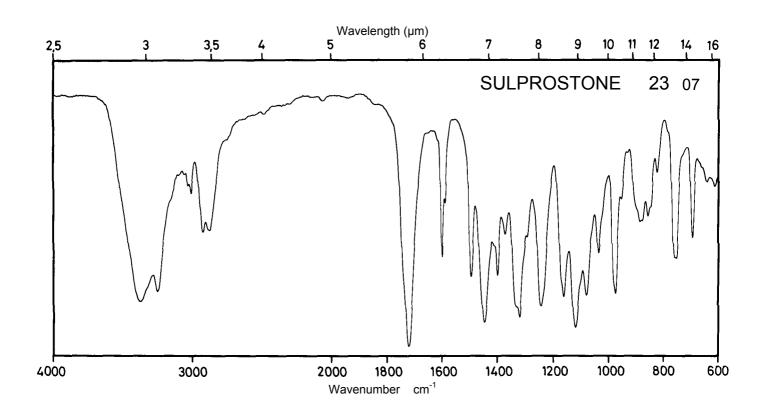

Name STANOZOLOL

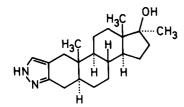



 M_{r}

328.5


Concentration 3 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	224 nm		230 nm	
E ^{1%} 1cm	146		203	
3	4800		6650	

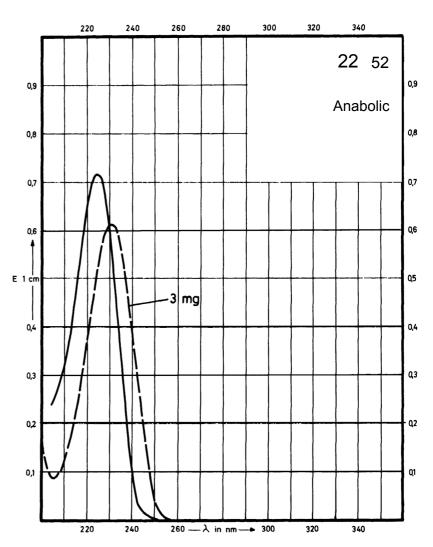

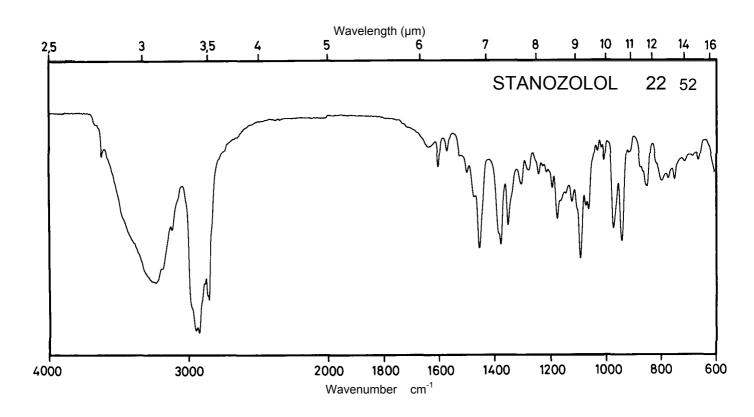





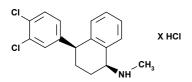


 M_{r} 342.7

Concentration 1 mg / 100 ml 25 mg / 100 ml

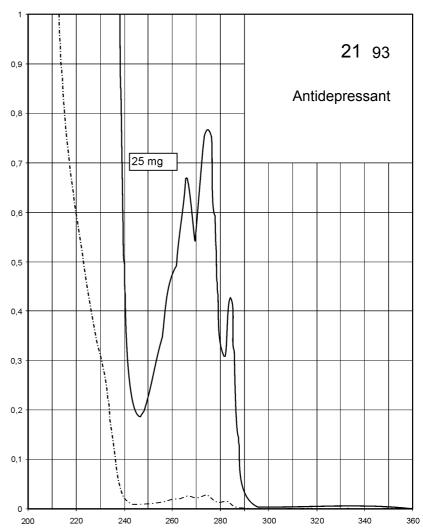

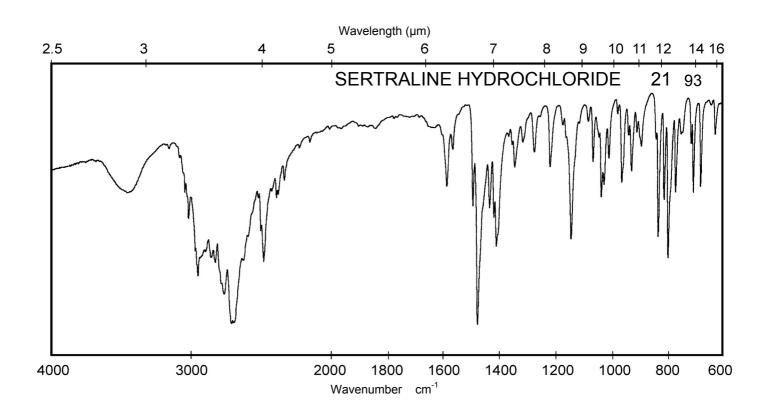




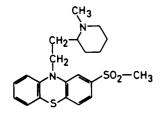


402.6


Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	315 nm 266 nm		311 nm 265 nm	266 nm
E ^{1%} 1cm	81 792		75 700	659
ε	3250 31890		3000 28180	26520

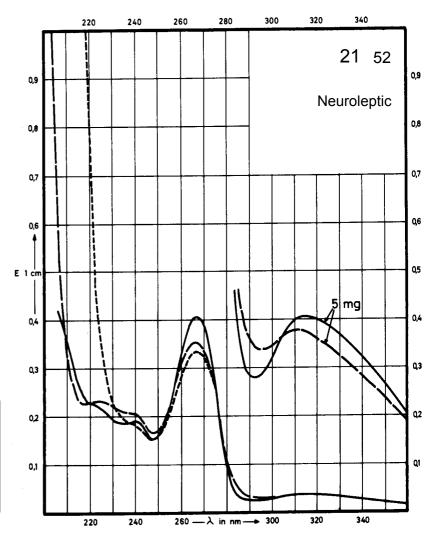

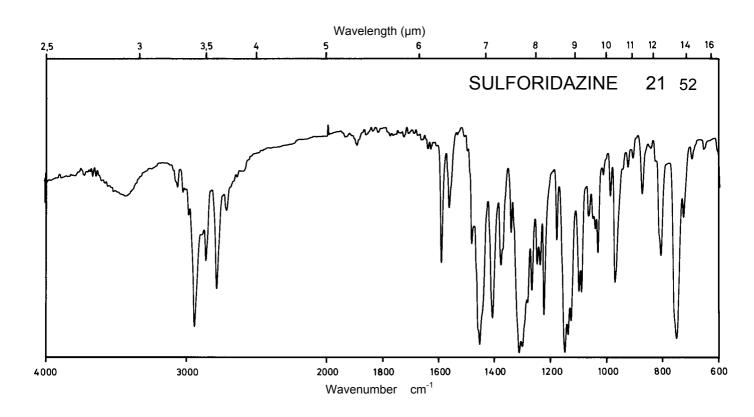


Name SULPIRIDE

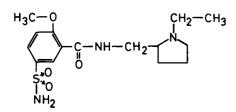


 M_{r}


341.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	289 nm		292 nm	292 nm
E ^{1%} 1cm	70		70	64
3	2390		2390	2180

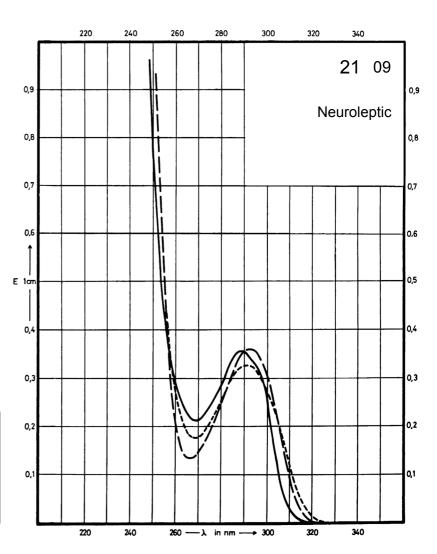

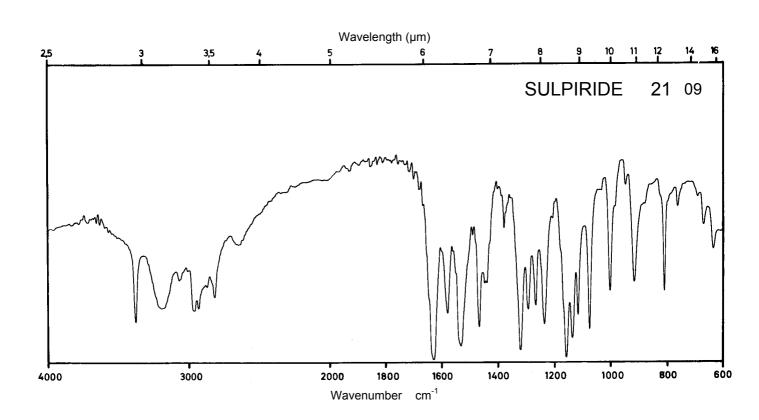




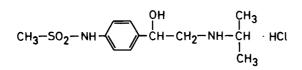


 M_r 308.8


Concentration 1 mg / 100 ml 20 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	274 nm 232 nm		269 nm	250 nm
E ^{1%} 1cm	18 472		13.7	487
3	555 14560		420	15020

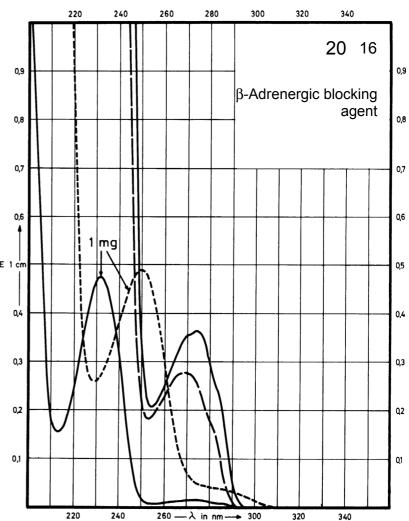

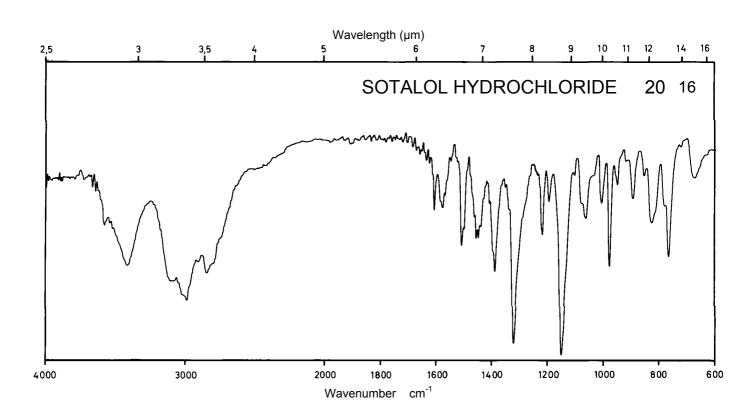


Name SCOPOLAMINE HYDROBROMIDE

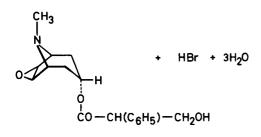



 M_{r}


438.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	257 nm		257 nm	
E ^{1%} 1cm	3.2		5.1	
3	140		220	

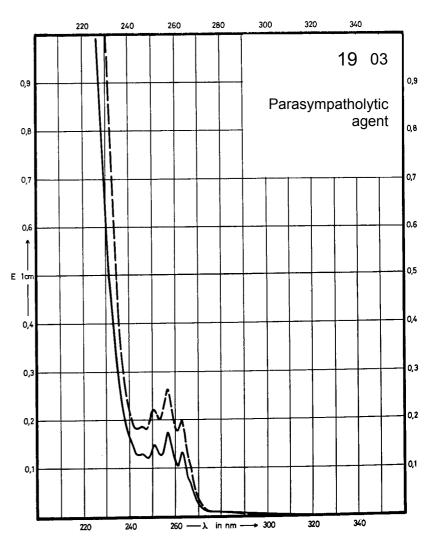


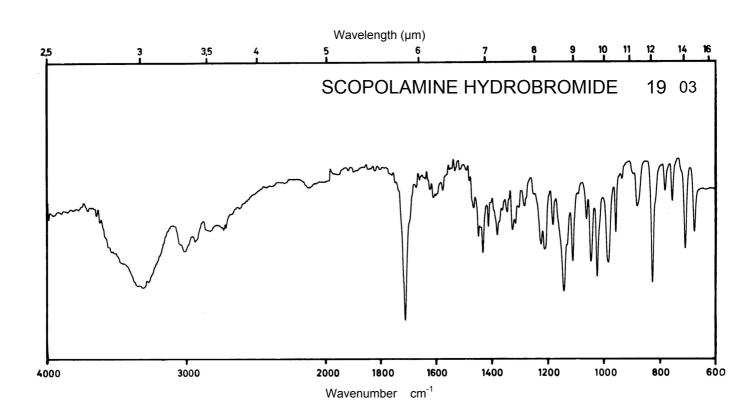

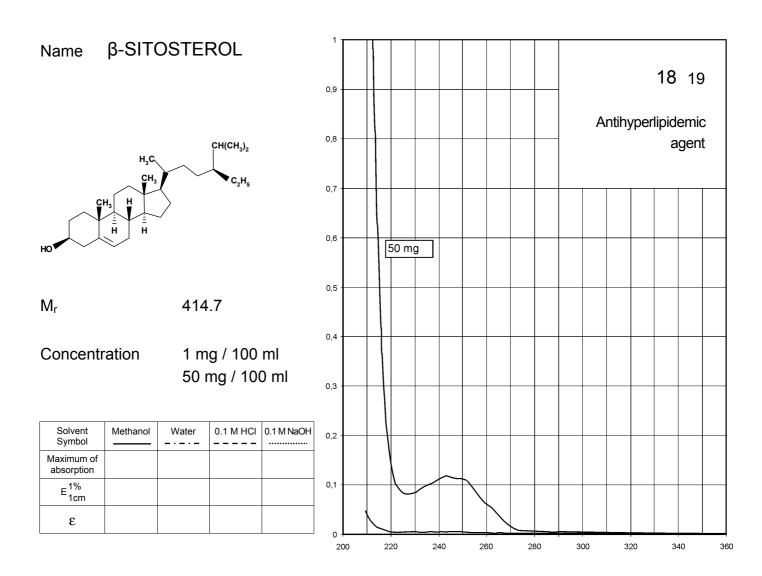

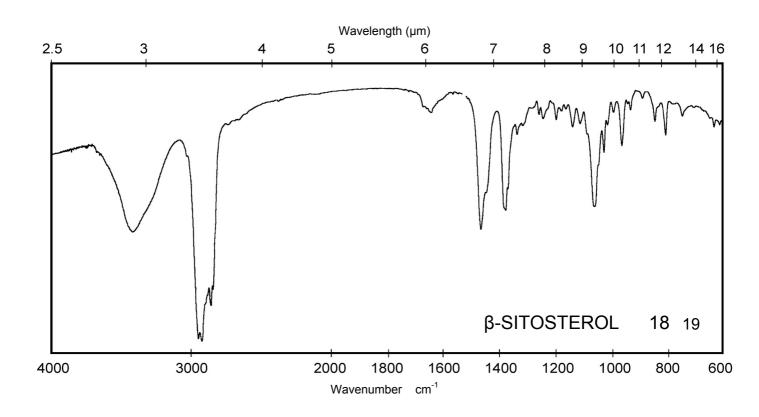





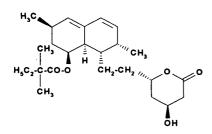


418.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	Decom- position observed	247 nm 239 nm 231 nm	247 nm 239 nm 231 nm	247 nm 239 nm 231 nm
E ^{1%} 1cm		332 498 448	350 524 474	364 551 514
ε		13900 20800 18800	14600 21900 19800	15200 23000 21500

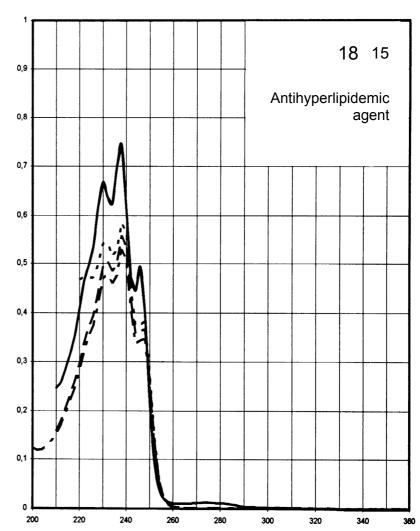

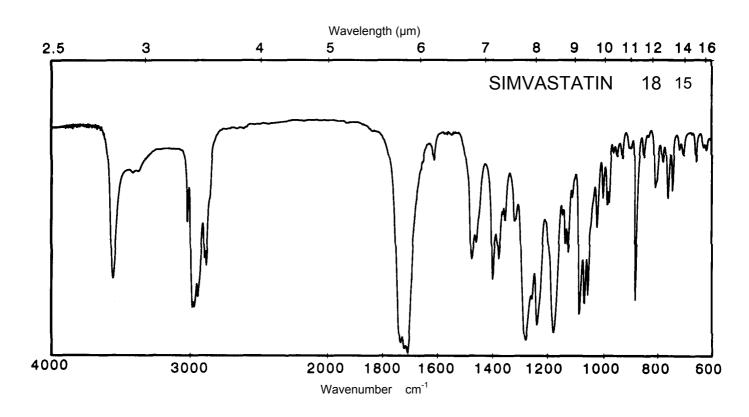





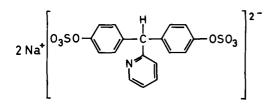


481.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	263 nm		264 nm	262 nm
E ^{1%} 1cm	123		169	119
3	5900		8120	5730

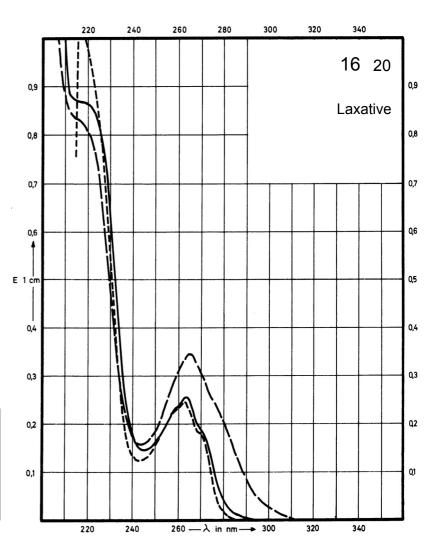

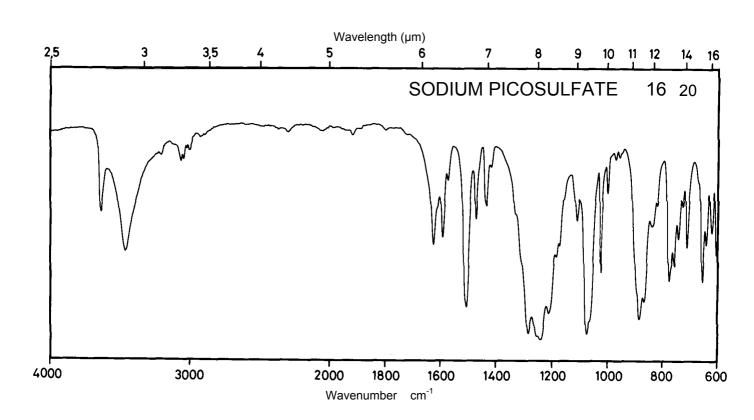



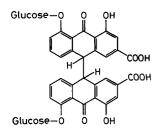


862.8


Concentration


5.2 mg / 100 ml


Solvent Symbol	Methanol	Water*	0.1 M HCI	0.1 M NaOH
Maximum of absorption		333 nm 270 nm		
E ^{1%} 1cm		150 200		
3		13000 17300		

^{* 1/15} M phosphate buffer pH 7.5

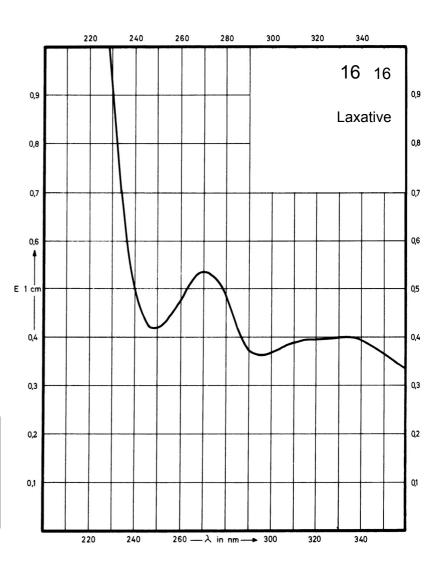


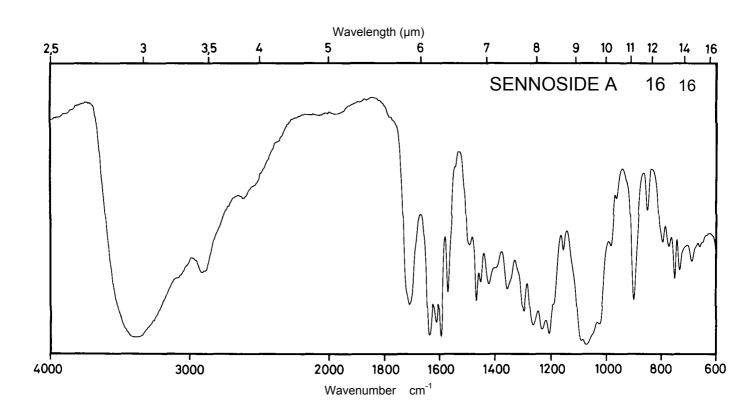

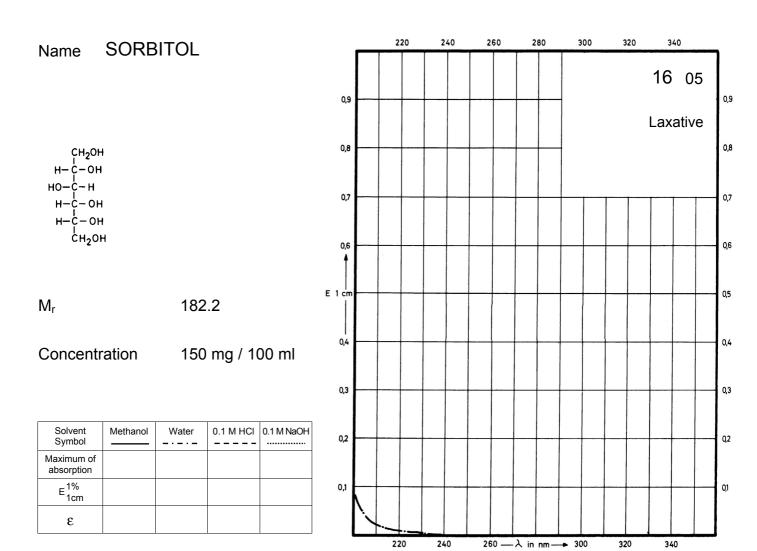

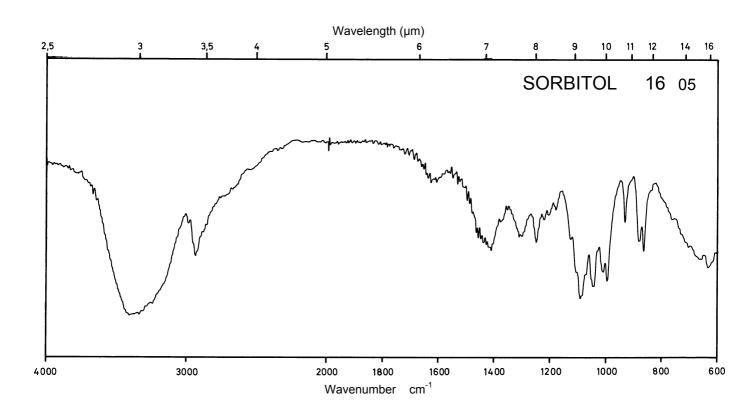





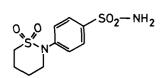


290.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	245 nm		244 nm	240 nm
E ^{1%} _{1cm}	393		343	372
3	11400		10000	10800

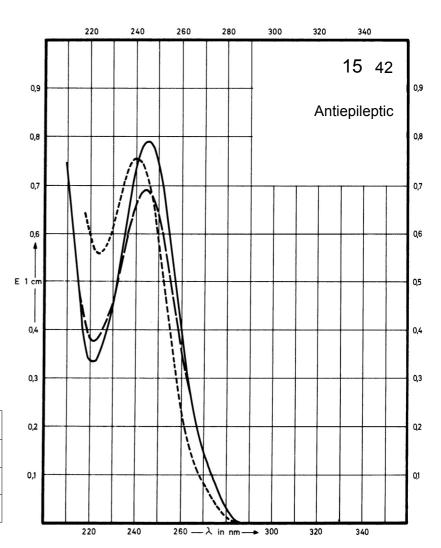

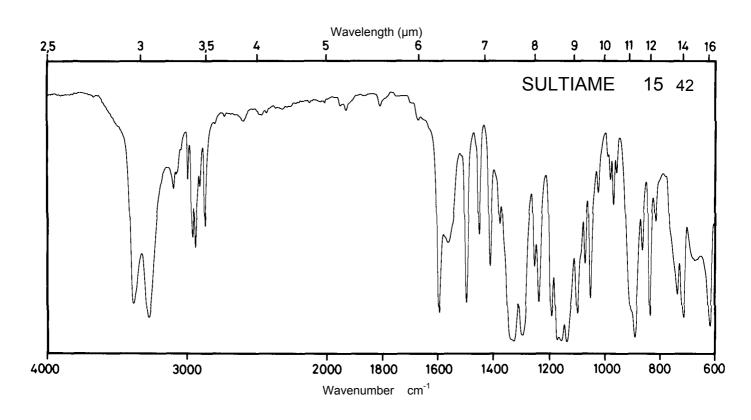





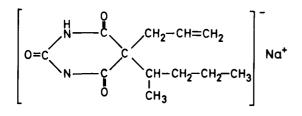


 M_r 260.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				243 nm
E ^{1%} 1cm				248
ε				6460

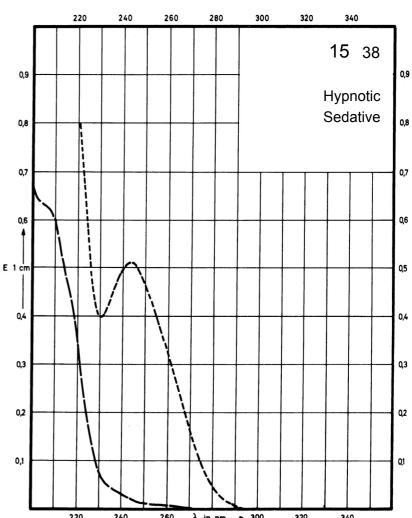

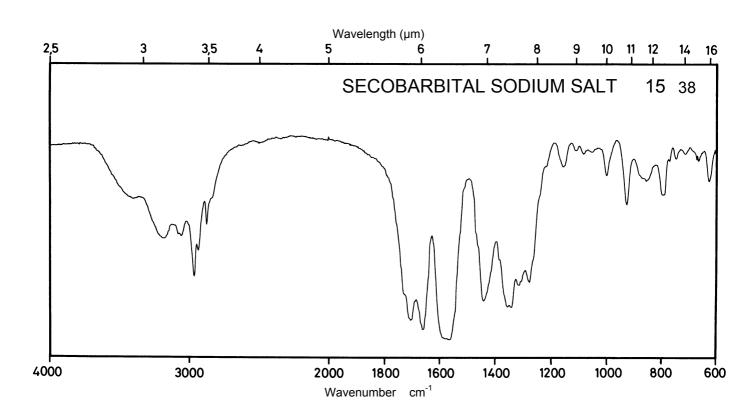





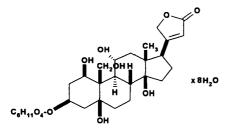


728.8

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	217 nm	220 nm	220 nm	
E ^{1%} 1cm	206	219	220	
3	15000	16000	16000	

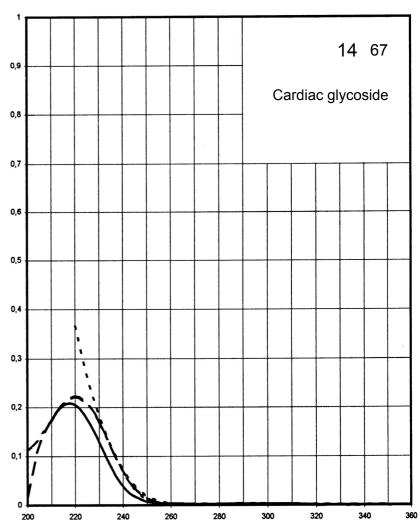

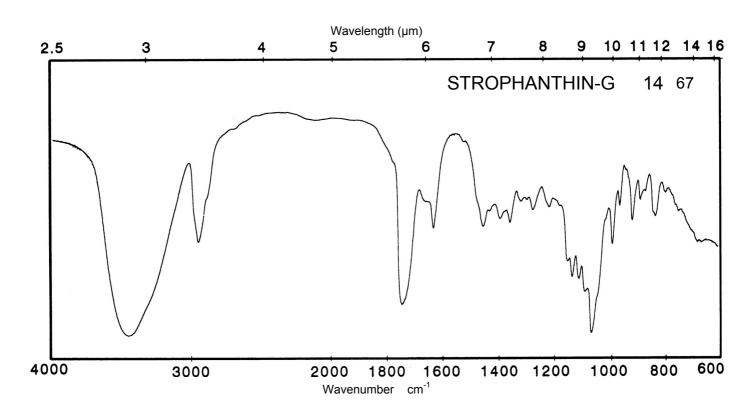





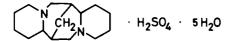


 M_{r} 422.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

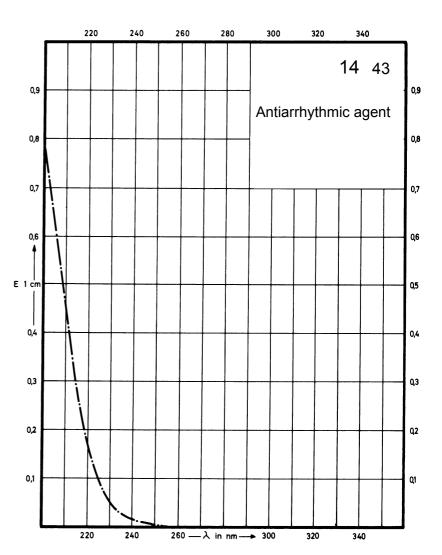

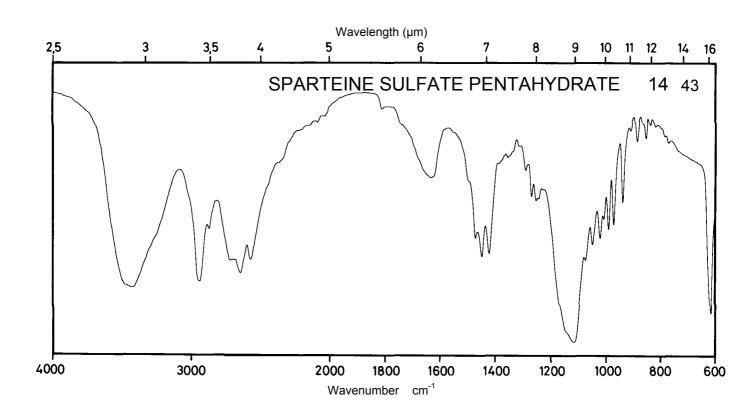





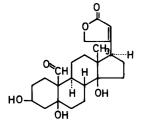

404.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
3				

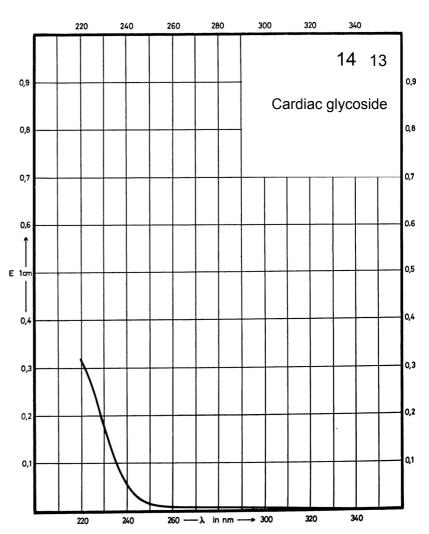

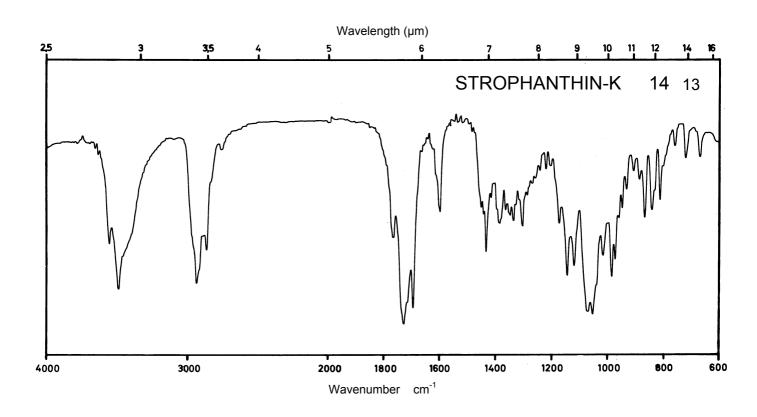



 $CH_3-CH=CH-CH=CH-COOK$


 M_{r}


150.2

Concentration

0.3 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	247 nm	254 nm	262 nm	254 nm
E ^{1%} _{1cm}	1710	1660	1750	1670
3	25700	24900	26300	25100

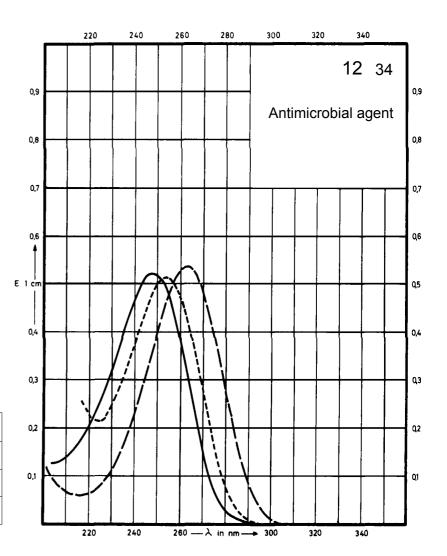

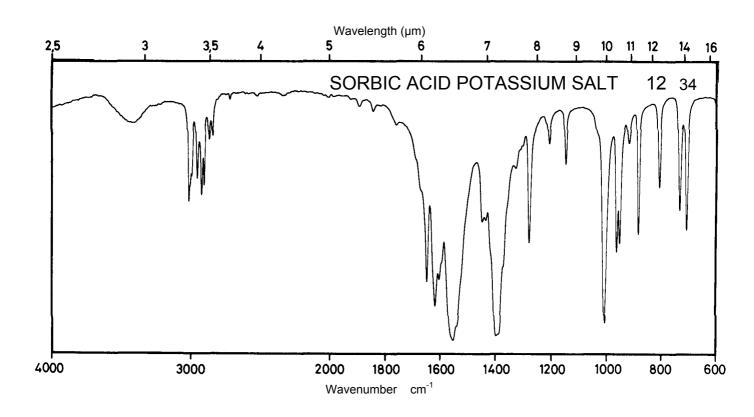



 $CH_3 - CH = CH - CH = CH - COOH$


 M_{r}


112.1

Concentration


0.2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	252 nm		262 nm	254 nm
E ^{1%} _{1cm}	2335		2306	2211
3	26180		25850	24790

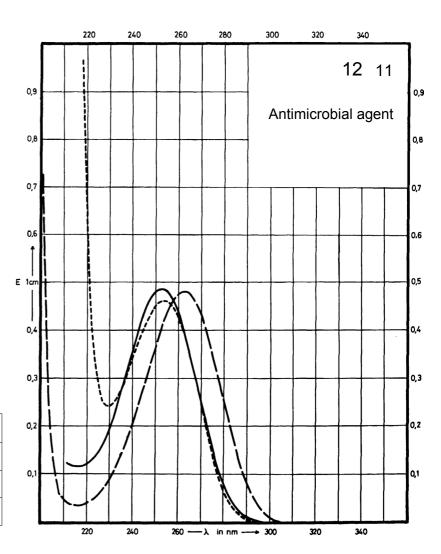

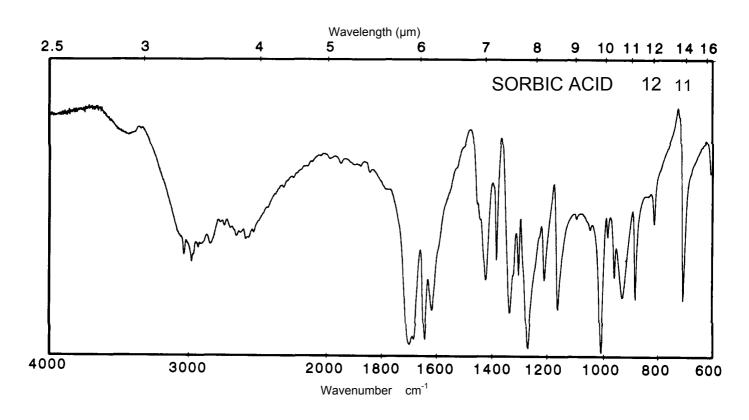



Name SODIUM TIMERFONATE

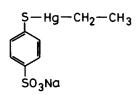

 M_r 440.9

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm	271 nm	255 nm	278 nm
E ^{1%} 1cm	300	323	325	315
ε	13210	14240	14330	13890

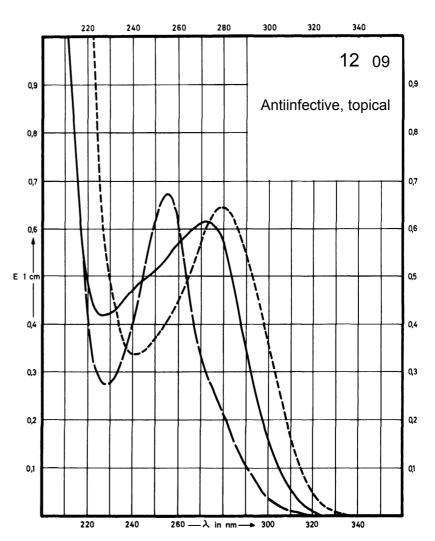

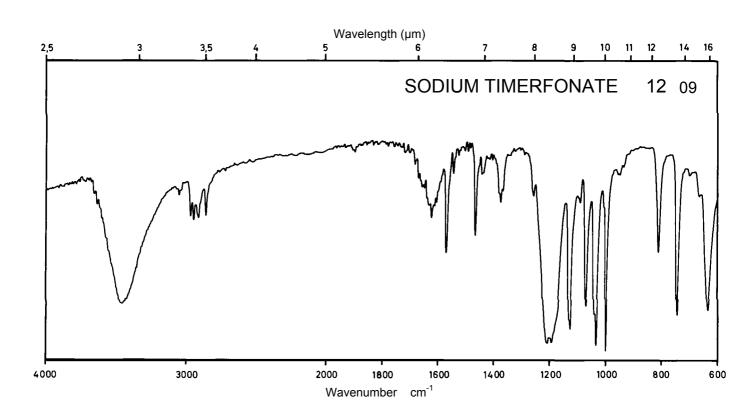


Name SALICYLIC ACID SODIUM SALT





 M_{r}


160.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	297 nm 228 nm	296 nm 230 nm	302 nm 236 nm	296 nm
E ^{1%} 1cm	242 410	225 423	229 536	220
3	3880 6600	3600 6770	3660 8600	3520

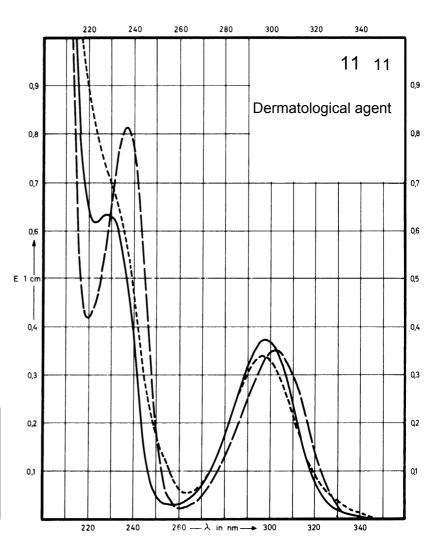

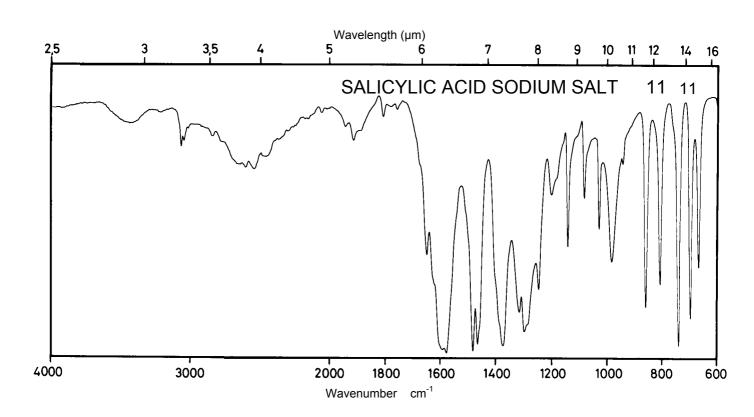

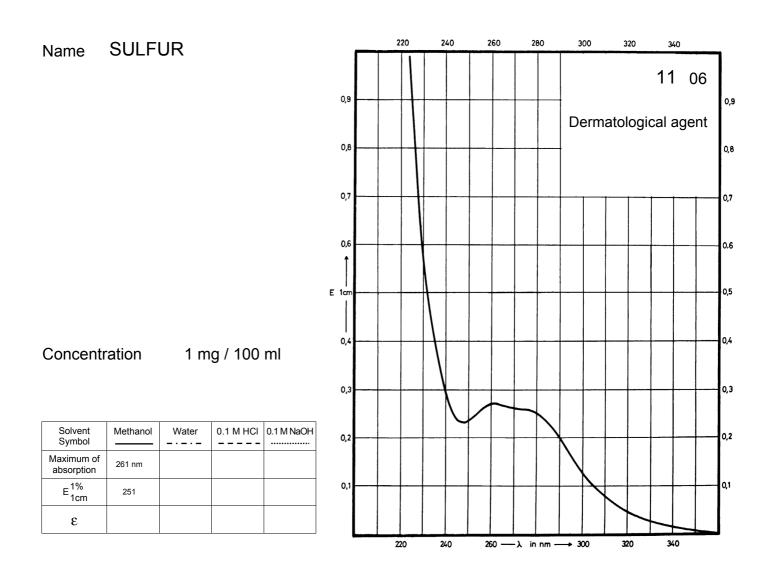

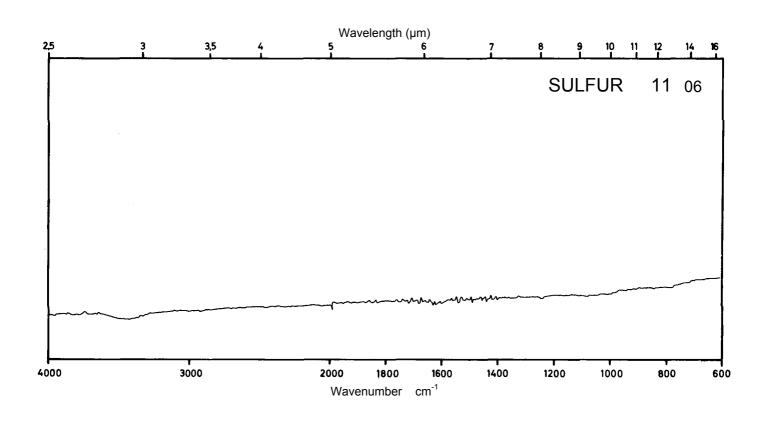


Name SALICYLIC ACID


СООН

 M_{r}


138.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	302 nm 234 nm		303 nm 237 nm	296 nm
E ^{1%} 1cm	285 547		262 613	254
3	3940 7550		3620 8470	3510

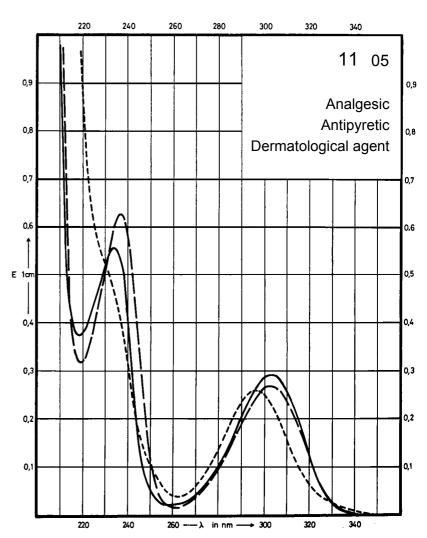

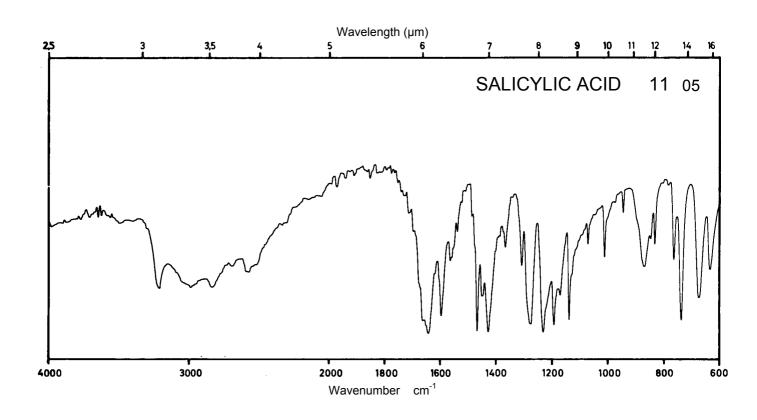

Name SULFAPHENAZOLE

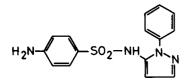



 M_{r}

314.4


Concentration 1 mg / 100 ml 2 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm			249 nm
E ^{1%} 1cm	652			700
ε	20500			22000

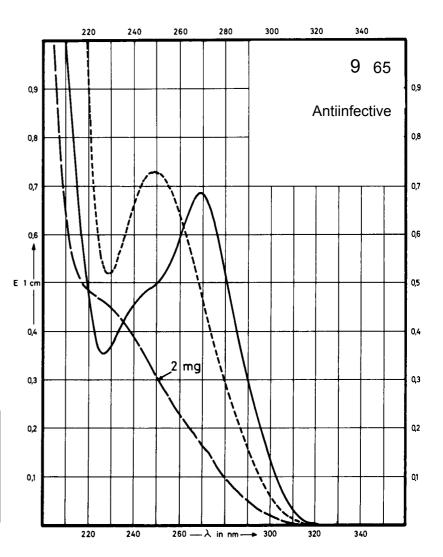

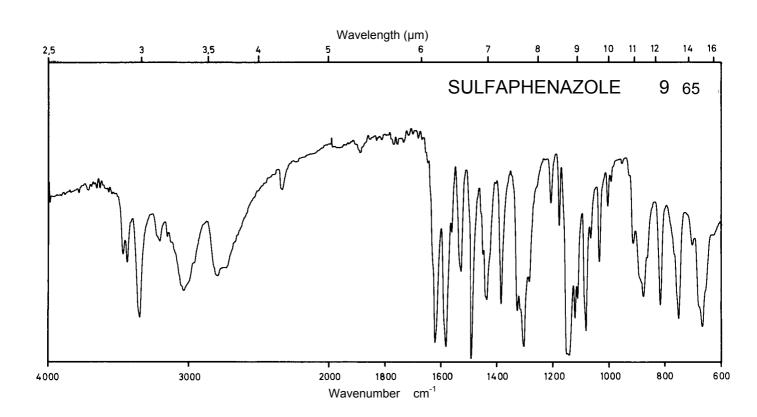





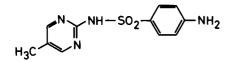


264.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		318 nm 244 nm	243 nm
E ^{1%} 1cm	800		115 542	900
3	21120		3040 14320	23800

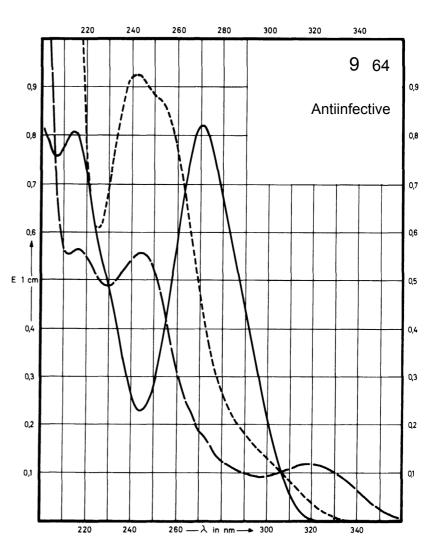

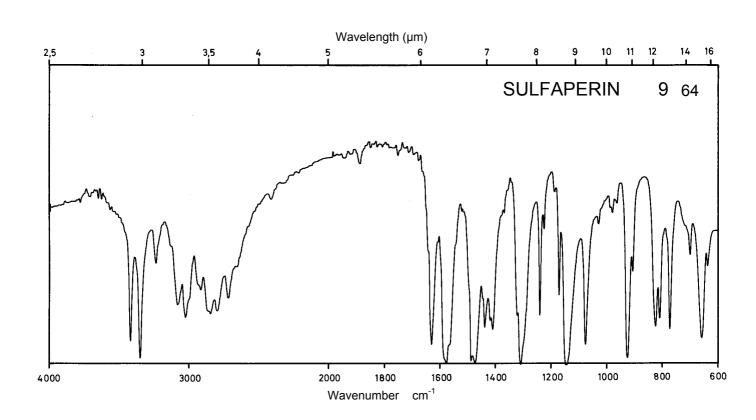





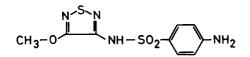


 M_r 286.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm	265 nm	288 nm	301 nm 257 nm
E ^{1%} 1cm	770	606	406	392 586
3	22040	17360	11640	11220 16770

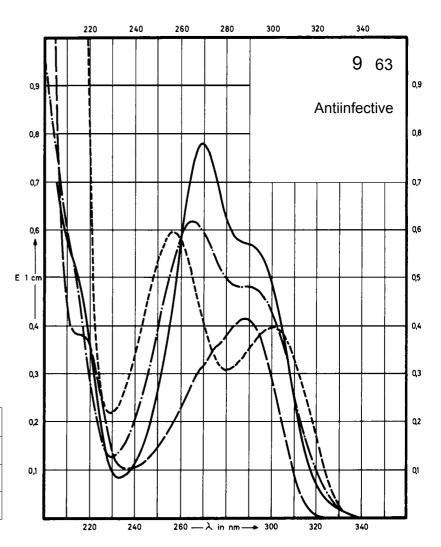

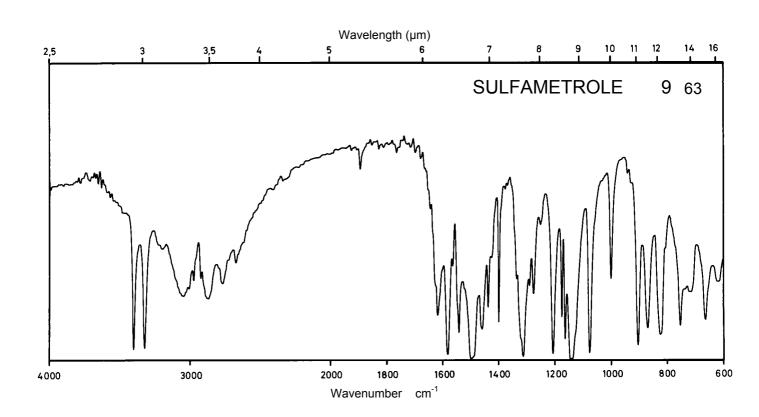




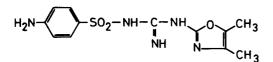


 M_{r}


309.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm		264 nm 224 nm	280 nm 259 nm
E ^{1%} 1cm	495		258 193	351 362
ε	15300		7970 5980	10850 11200

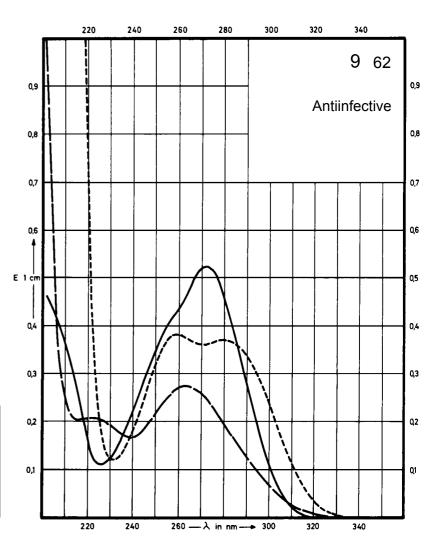

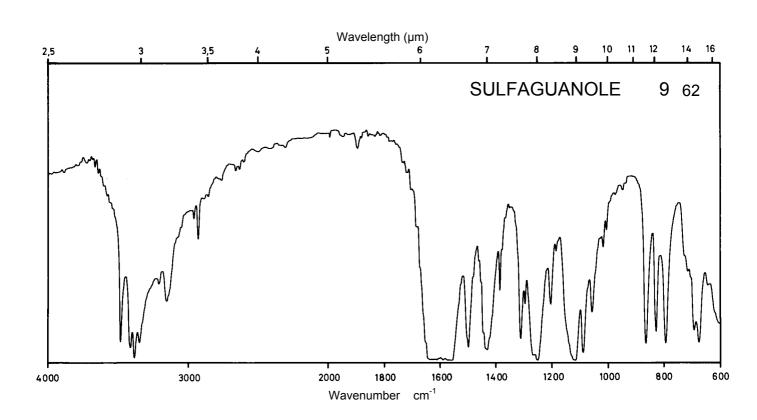




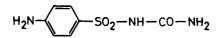


 M_r 215.2


Concentration 0.5 mg / 100 ml 1 mg / 100 ml

Solvent Methanol Water 0.1 M HCI 0.1 M NaOH Symbol Maximum of absorption 266 nm 218 nm 270 nm 255 nm 200 483 924 773 4300 10400 16640 ε 19880

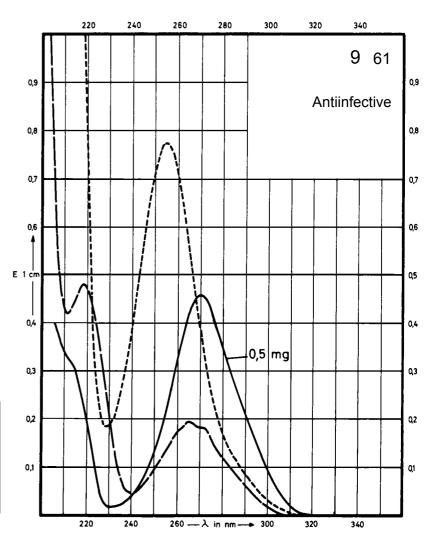

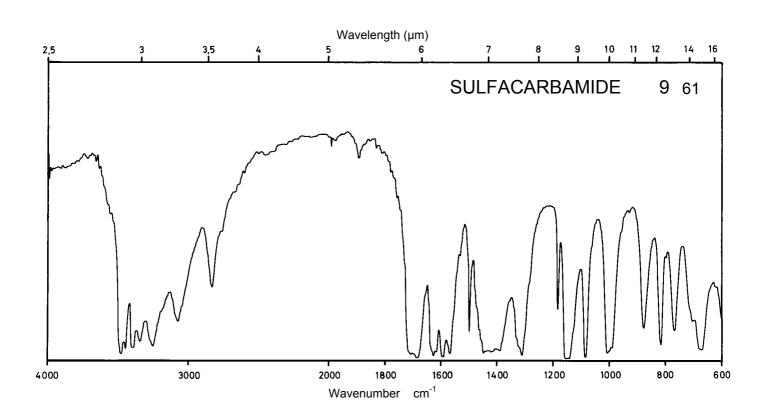




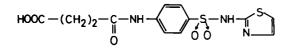


 M_r


355.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	285 nm 258 nm		281 nm 258 nm	257 nm
E ^{1%} 1cm	640 629		590 593	660
ε	22750 22350		20970 21080	23460

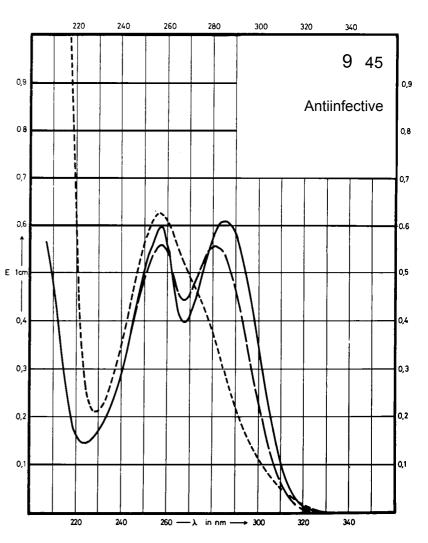

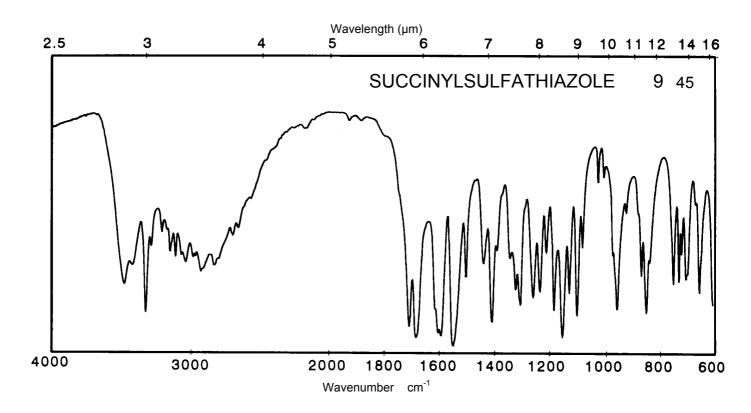


Name SULFAMETHOXAZOLE

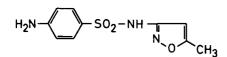


 M_{r}


253.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm		265 nm	256 nm
E ^{1%} 1cm	785		172	675
3	19880		4360	17100

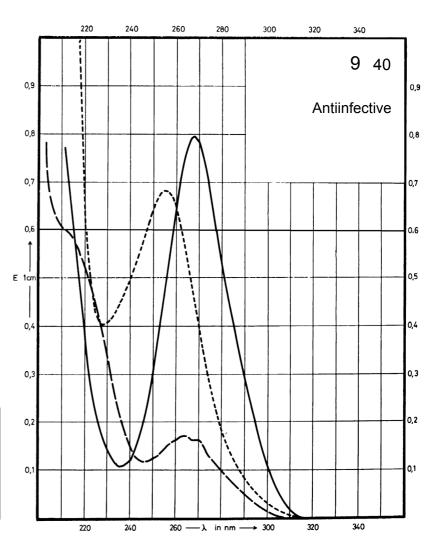

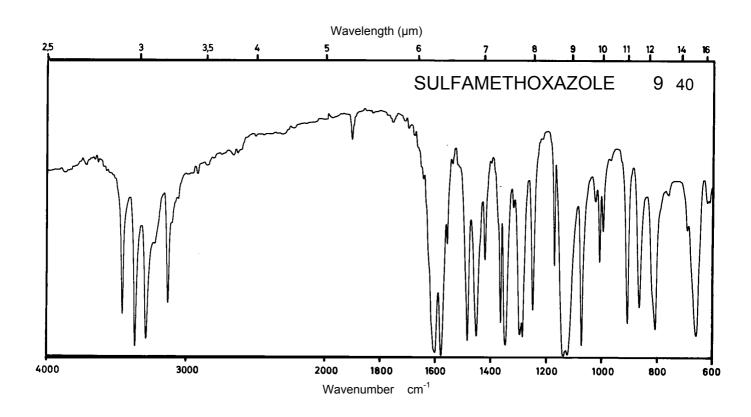





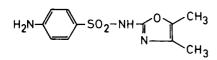

267.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm		249 nm	250 nm
E ^{1%} 1cm	840		375	685
ε	22420		10030	18310

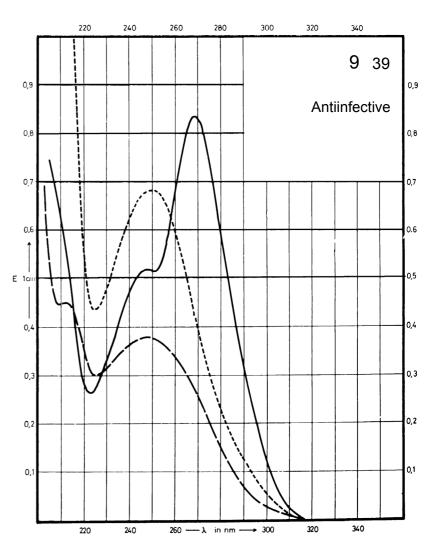

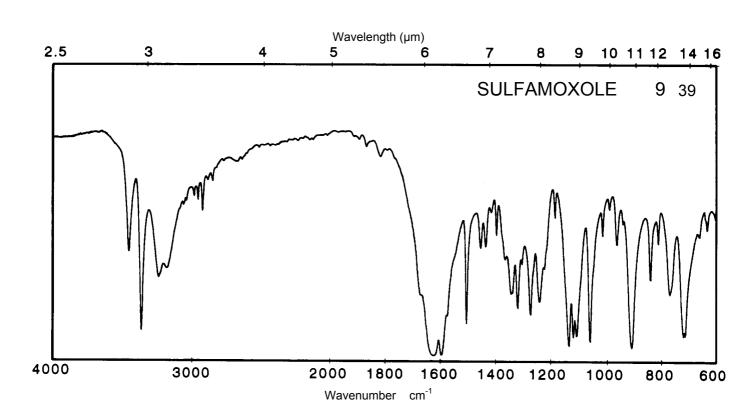


Name SULFISOMIDINE


 M_r

Concentration 1 mg / 100 ml


278.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm		262 nm	262 nm
E ^{1%} 1cm	729		648	764
3	20290		18030	21260

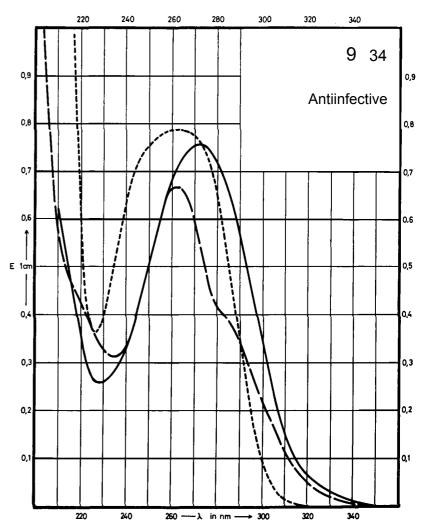

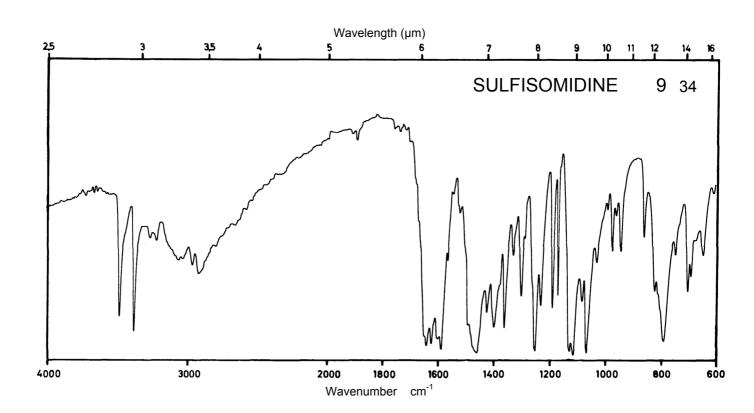



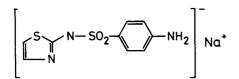



277.3

Concentration


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	287 nm 260 nm		280 nm	256 nm
E ^{1%} _{1cm}	731 603		462	661
ε	20300 16700		12800	18300

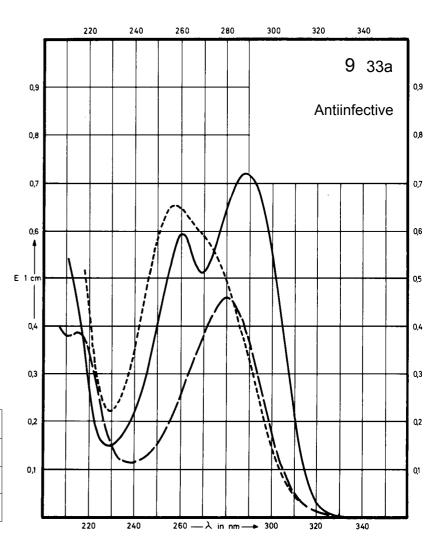

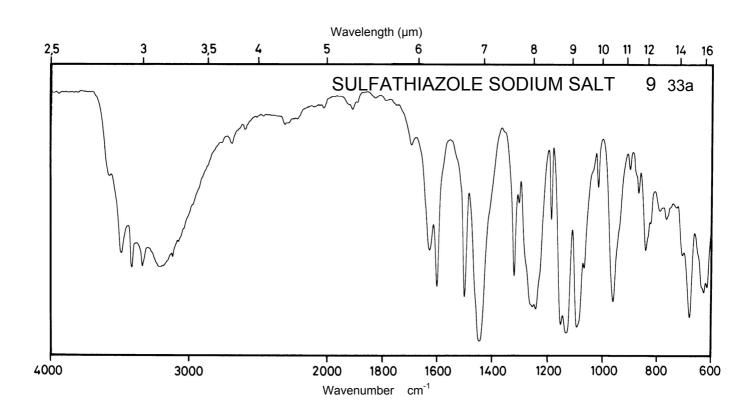

Name SULFATHIAZOLE

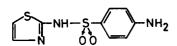


 M_{r}


255.3

Concentration


0.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	288 nm 258 nm		280 nm	256 nm
E ^{1%} 1cm	835 633		505	716
ε	21320 16160		12890	18280

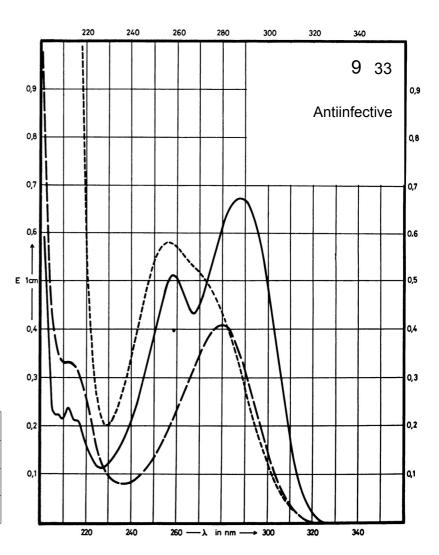

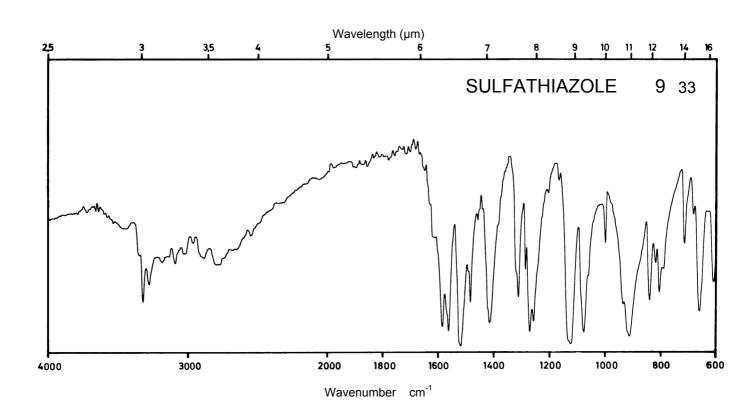



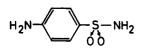



 M_{r}

172.2

Concentration 0.5 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	262 nm		269 nm 262 nm	250 nm
E ^{1%} 1cm	1116		88 112	932
ε	19220		1520 1930	16050

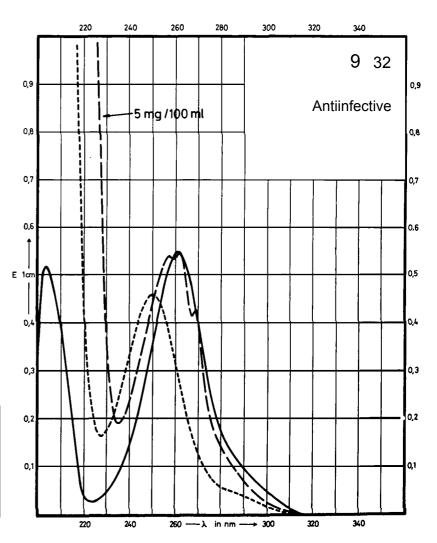

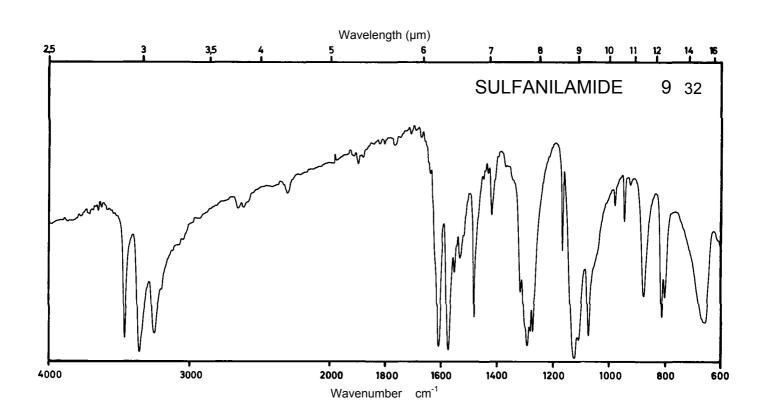

Name SULFAMETHOXY-PYRIDAZINE


 M_{r}


280.3

Concentration 1


1.1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm		316 nm	252 nm
E ^{1%} 1cm	775		163	723
ε	21720		4570	20270

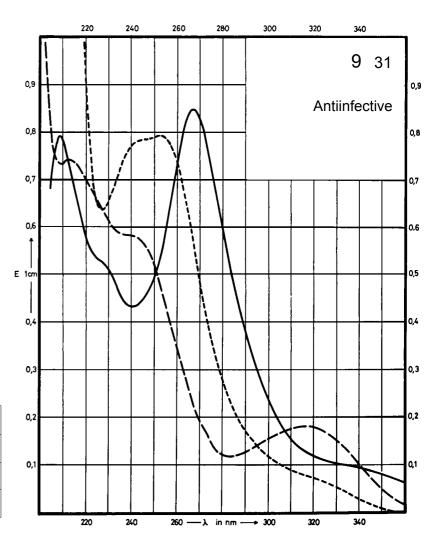

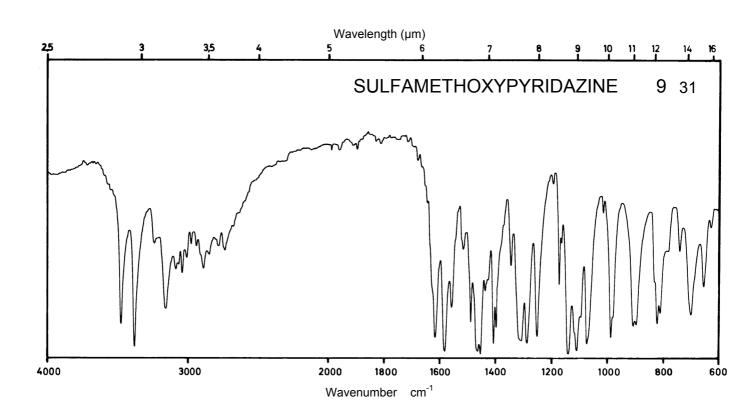


Name SULFAMETOXY-DIAZINE

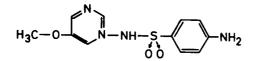


 M_{r}


280.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm 231 nm		218 nm	318 nm 244 nm
E 1% 1cm	739 573		608	105 932
ε	20710 16060		17040	2940 26120

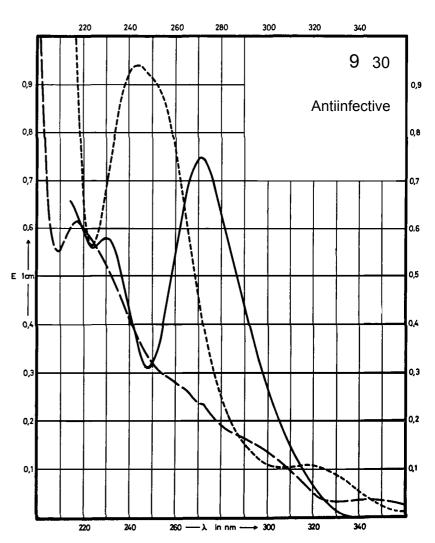

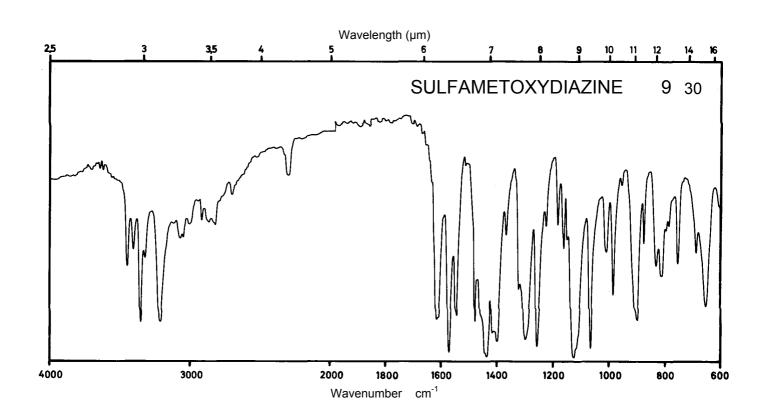

Name SULFAMERAZINE

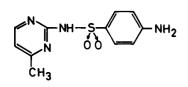


 M_{r}


264.3

Concentration

1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	270 nm		304 nm 242 nm	255 nm 241 nm
E 1% 1cm	860		198 596	828 820
ε	22730		5230 15750	21880 21670

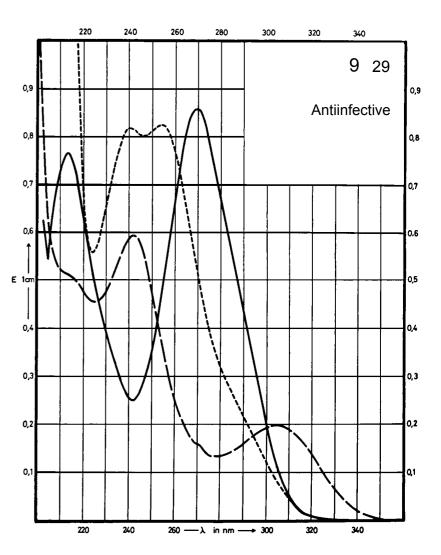

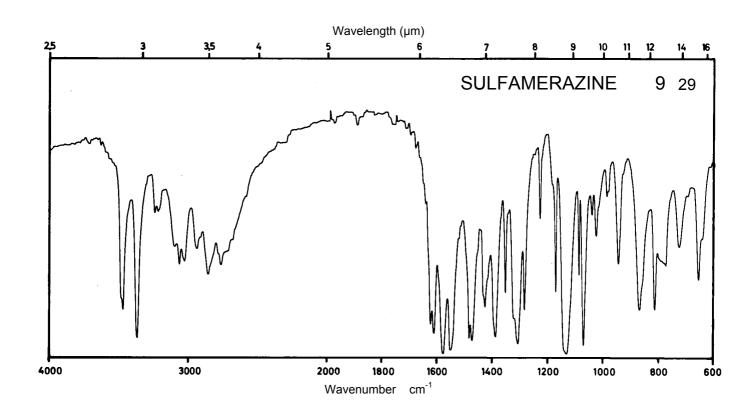



Name SULFAGUANIDINE

 M_{r}

214.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	264 nm		271 nm 264 nm	258 nm
E ^{1%} 1cm	879		95 103	738
ε	18830		2030 2210	15810

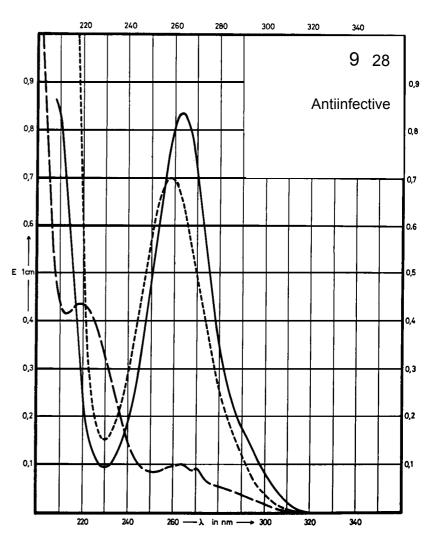

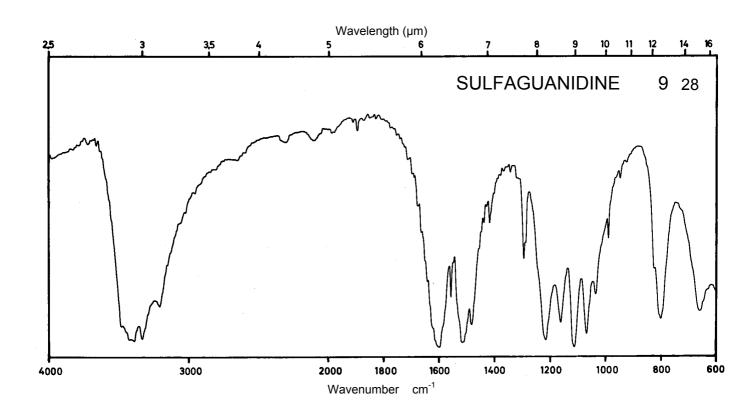



Name SULFAFURAZOLE

 M_{r}


267.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		265 nm	253 nm
E ^{1%} 1cm	720		235	796
ε	19250		6280	21280

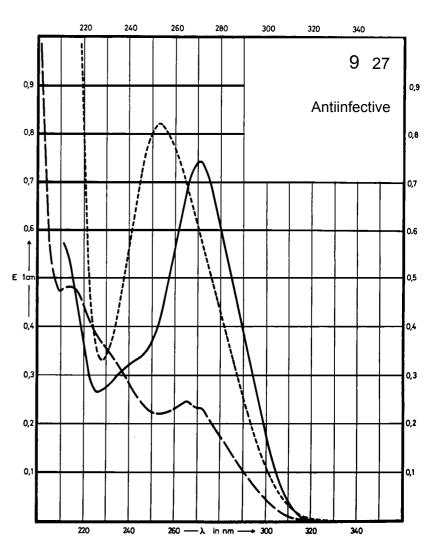

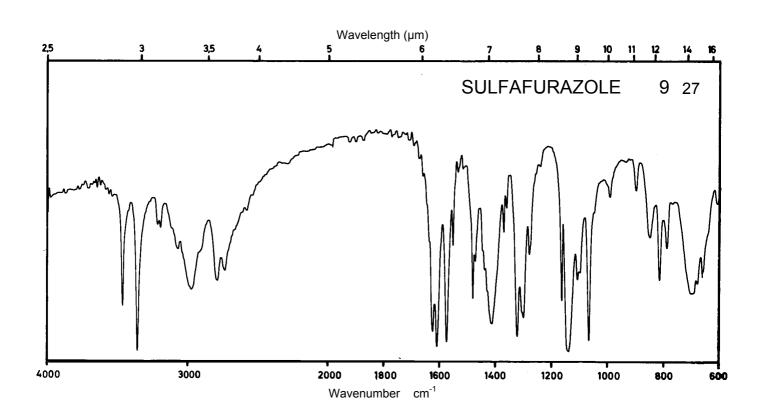





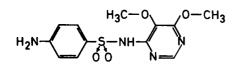


310.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	273 nm		264 nm	271 nm
E ^{1%} 1cm	744		406	775
3	23090		12600	24050

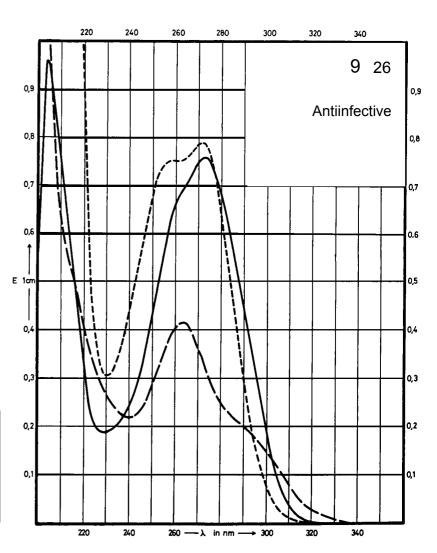

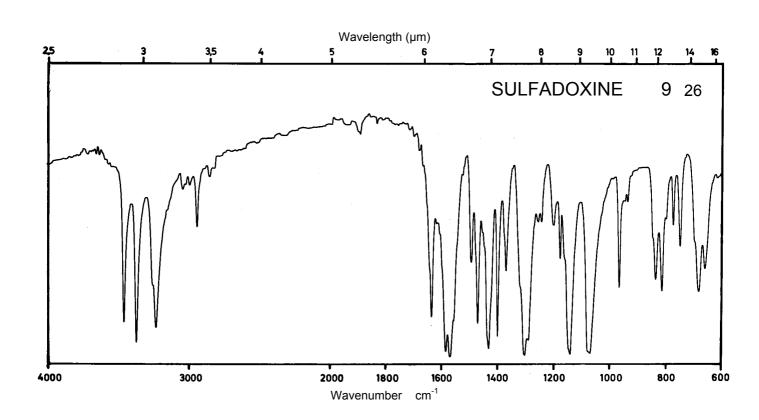



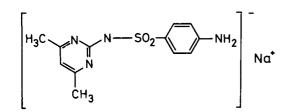



300.3

Concentration

0.95 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm		301 nm 243 nm	257 nm 243 nm
E ^{1%} 1cm	757		235 499	725 701
ε	22700		7060 15000	21800 21000

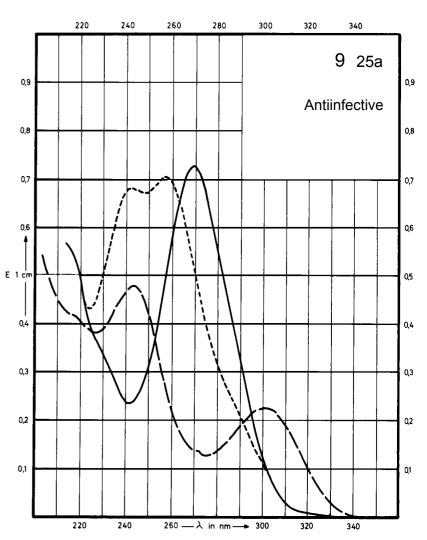

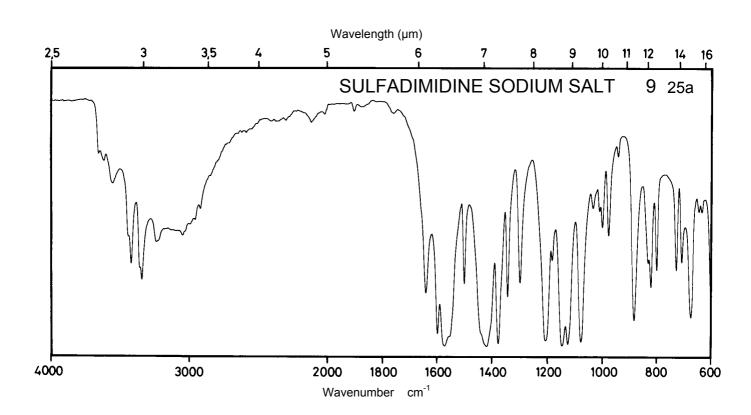

Name SULFADIMIDINE


 M_{r}


278.3

Concentration 1 mg


1 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	269 nm		301 nm 243 nm	257 nm 242 nm
E ^{1%} 1cm	812		248 531	780 756
ε	22600		6900 14780	21710 21040

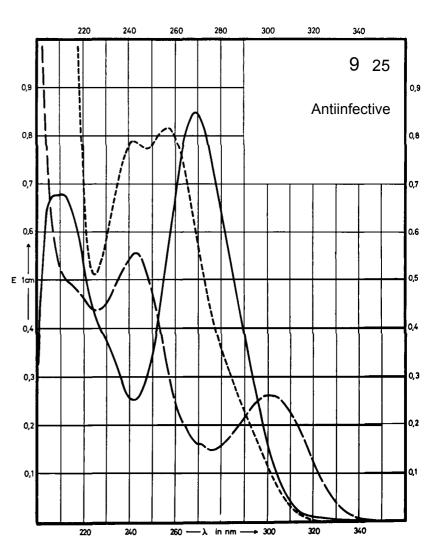

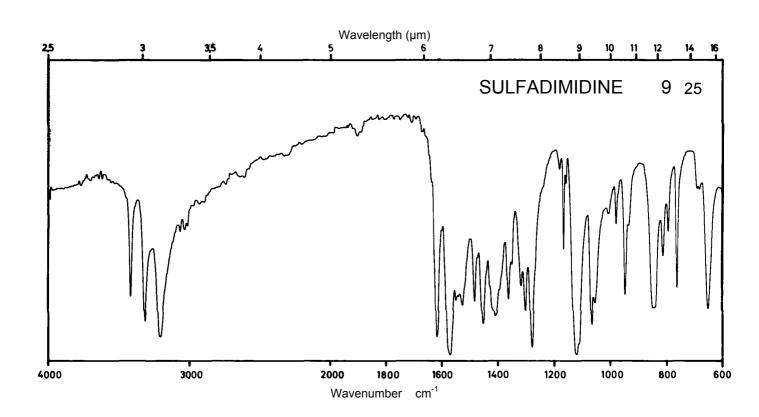




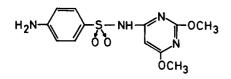


 M_r


310.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	272 nm		275 nm	269 nm
E ^{1%} 1cm	699		441	865
ε	21690		13680	26840

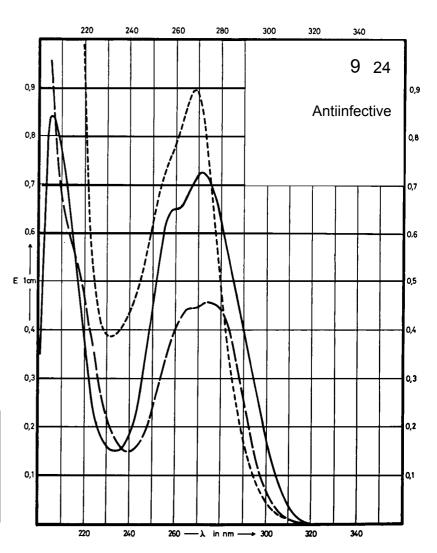

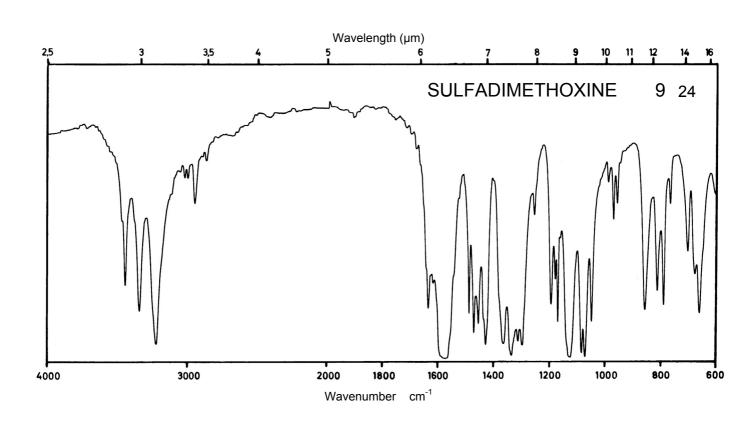

Name SULFADIAZINE

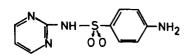


 M_{r}


250.3

Concentration


0.9 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	268 nm		305 nm 240 nm	254 nm 240 nm
E ^{1%} 1cm	851		116 511	861 869
ε	21300		2900 12790	21550 21750

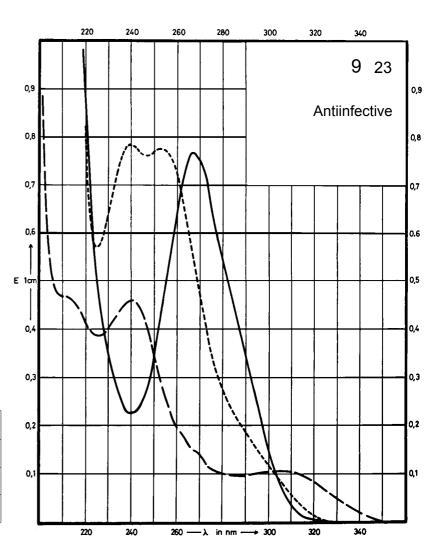

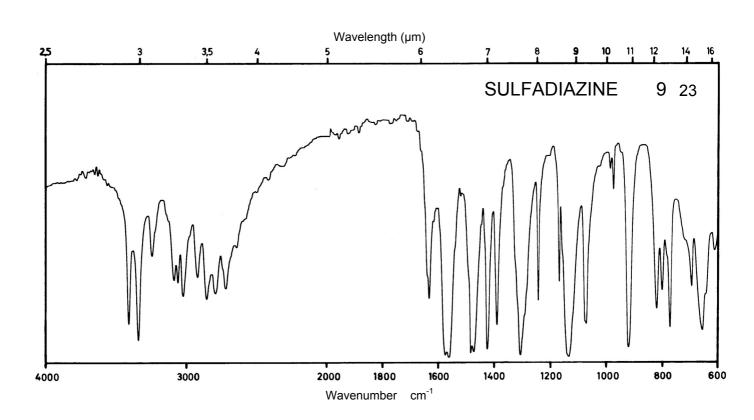





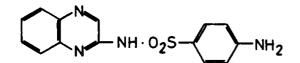

300.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	338 nm 267 nm 248 nm		348 nm 263 nm 248 nm	357 nm 250 nm
E ^{1%} 1cm	217 697 949		263 540 648	265 1069
3	6520 20930 28500		7900 16220 19460	7960 32100

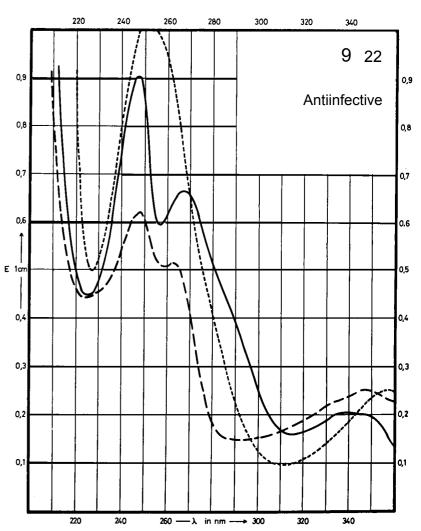

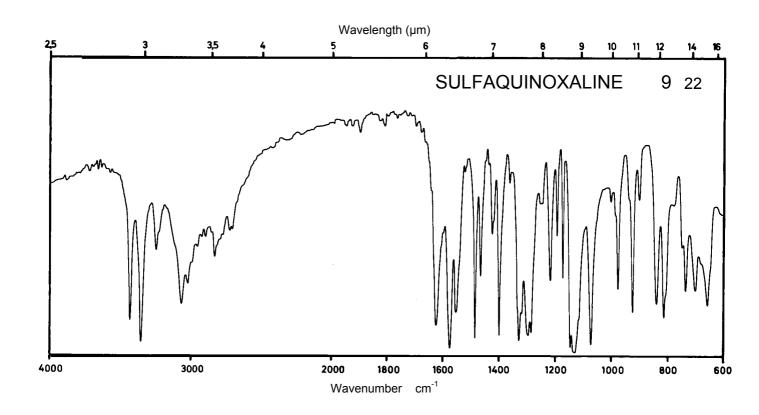


Name SULFACETAMIDE SODIUM SALT


$$\begin{bmatrix} SO_2 - N - CO - CH_3 \\ NH_2 \end{bmatrix} - Na^+$$

 M_r 254.2

Concentration 0.8 mg / 100 ml 1.5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	265 nm		271 nm 265 nm	256 nm
E 1% 1cm	683		197 194	650
ε	17400		5000 4950	16500

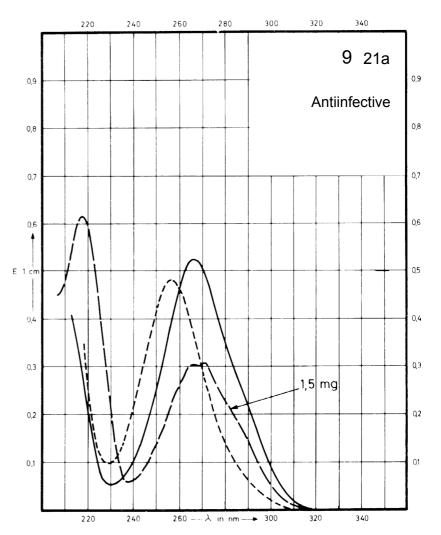

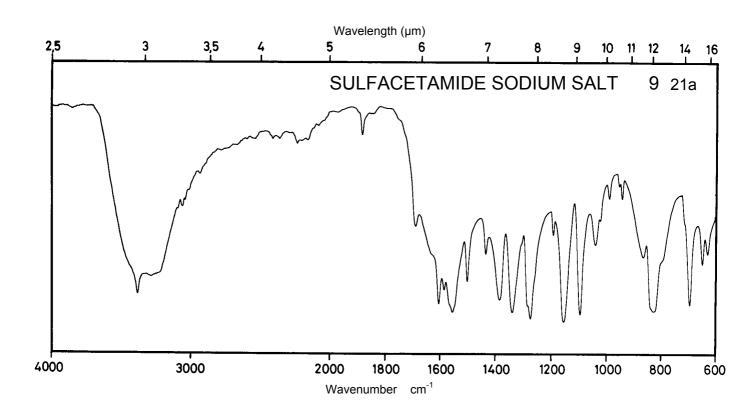

Name SULFACETAMIDE

SO₂-NH-CO-CH₃


 M_{r}


214.2

Concentration


0.8 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	271 nm		271 nm 265 nm	256 nm
E ^{1%} 1cm	920		246 244	650
ε	19700		5250 5200	16500

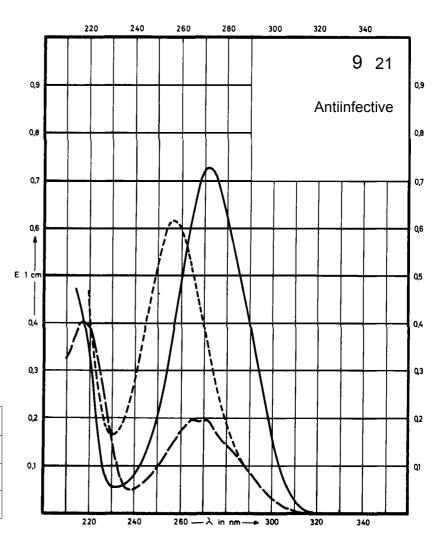

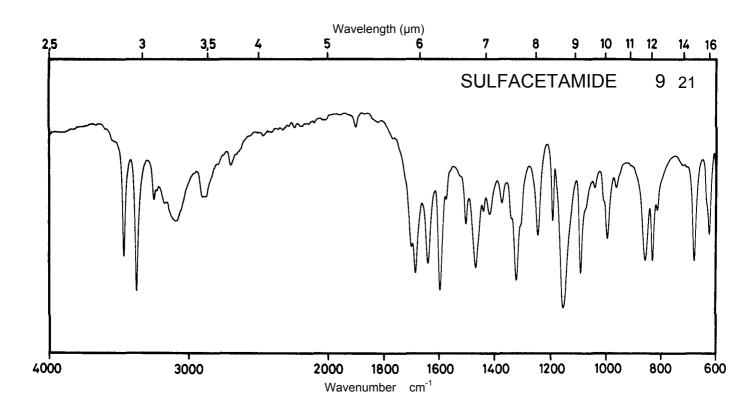



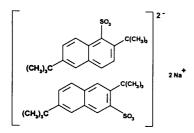



342.4

Concentration


0.5 mg / 100 ml

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of	323 nm	323 nm	323 nm	323 nm
	281 nm	283 nm	281 nm	281 nm
absorption	231 nm	232 nm	231 nm	231 nm
1%	31	39	39	38
E ^{1%} 1cm	158	149	149	146
1cm	2170	2000	2000	1812
	1070	1335	1300	1300
3	5400	5100	5100	5000
	74300	65000	65000	62100

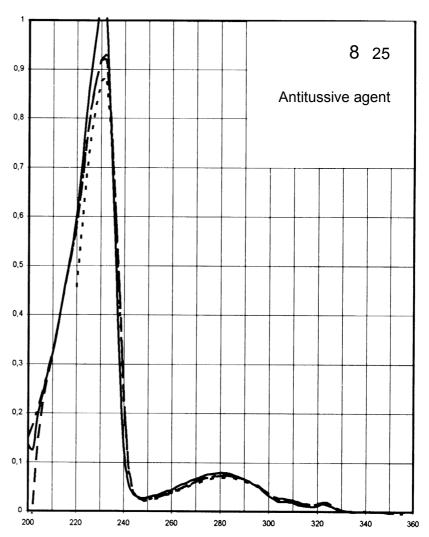

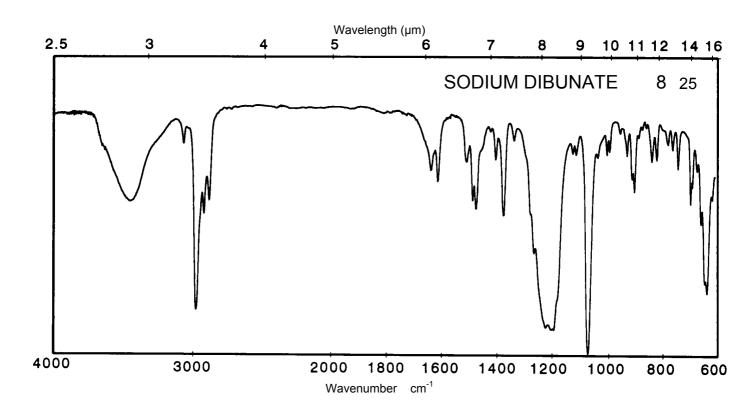





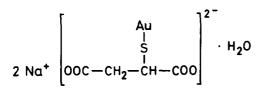


408.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} 1cm				
ε				

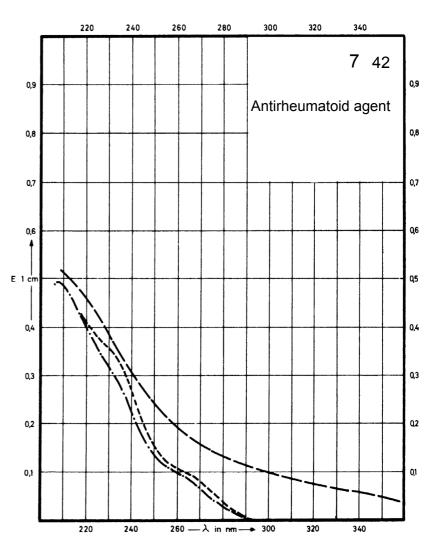

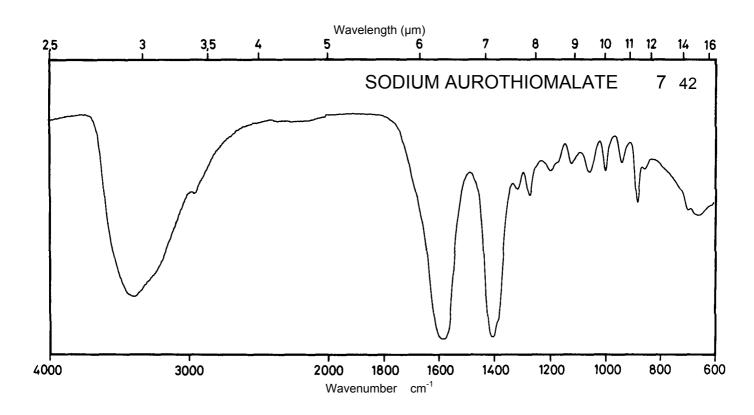





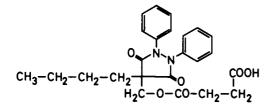

 M_r 438.5

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	237 nm		235 nm	263 nm
E ^{1%} 1cm	375		335	438
ε	16450		14700	19210

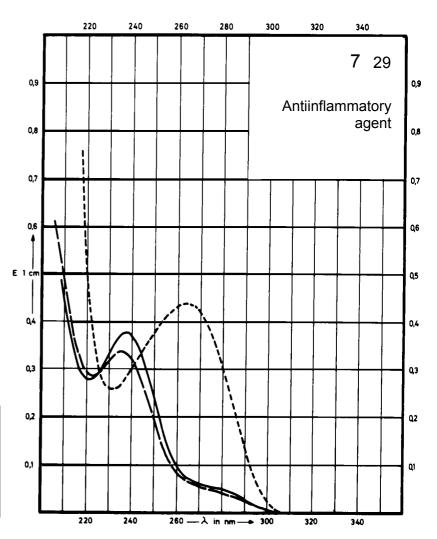


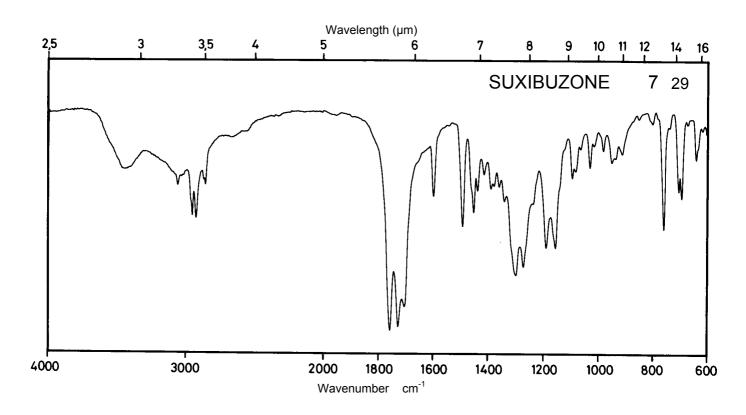


Name SALICYLIC ACID
TRICHLORISOBUTYL
ESTER


 M_{r} 297.6

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	308 nm 240 nm			
E ^{1%} 1cm	155 345			
ε	4600 10300			

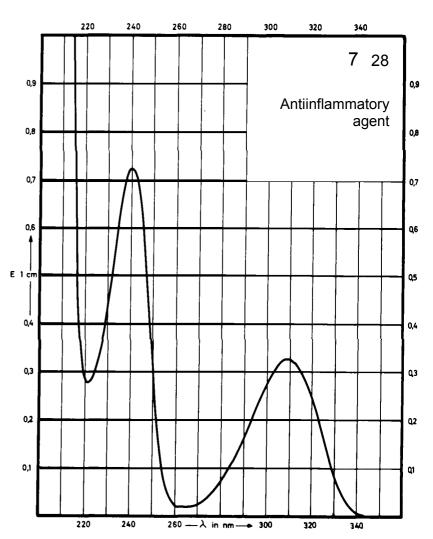

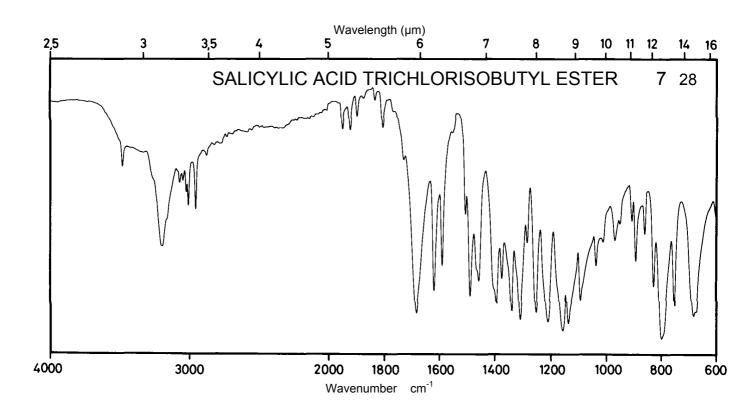





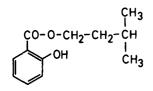

208.3

Solvent	Methanol	Water	0.1 M HCI	0.1 M NaOH
Symbol				
Maximum of absorption	304 nm 236 nm			
E ^{1%} 1cm	215 455			
ε	4480 9480			

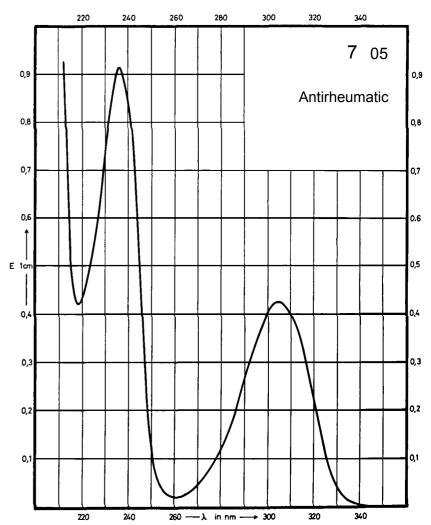

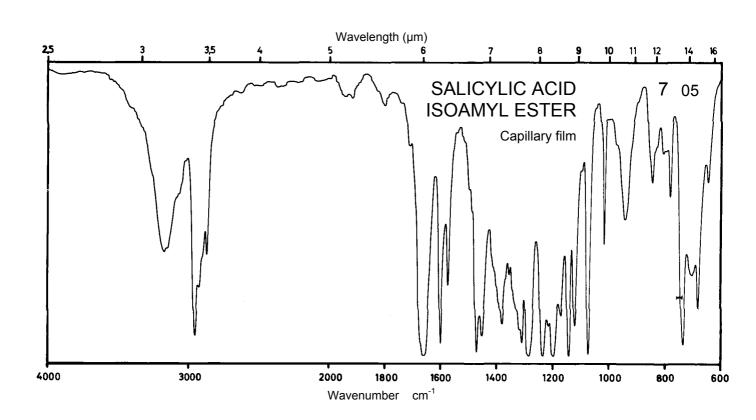





 $Na_2[Fe(CN)_5NO] \cdot 2 H_2O$

 M_r 298.0


Concentration 1 mg / 100 ml 5 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption		397 nm		Decom- position observed
E ^{1%} 1cm		0.70		
ε		21		

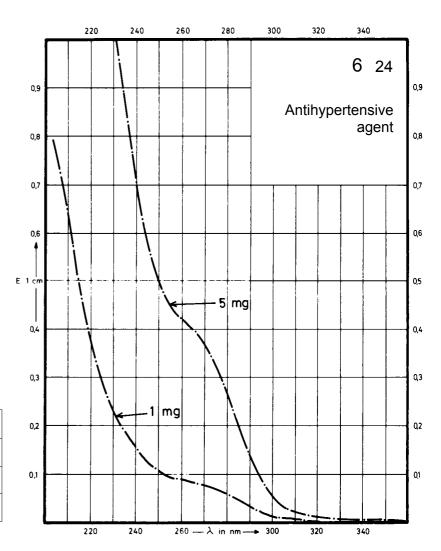

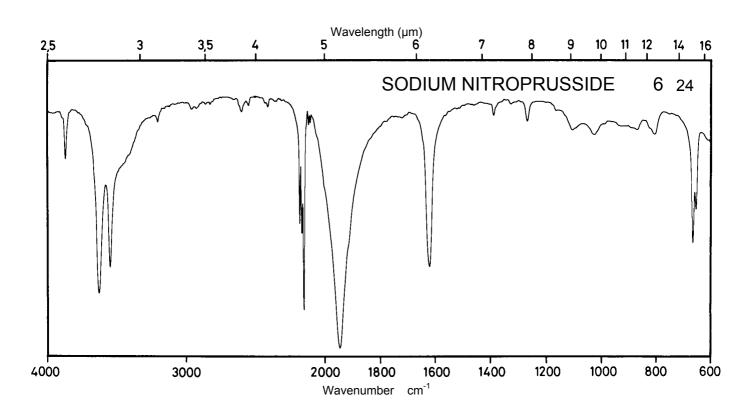





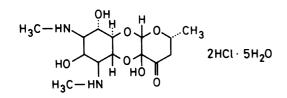



 M_{r} 495.4


Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption				
E ^{1%} _{1cm}				
ε				

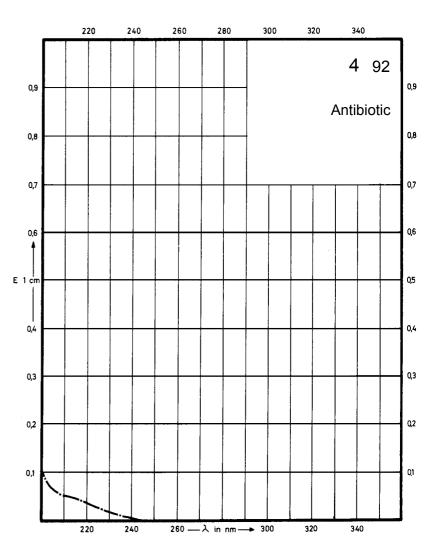


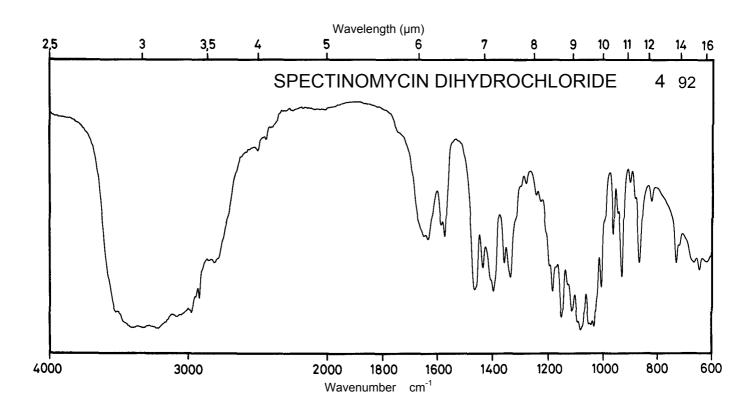

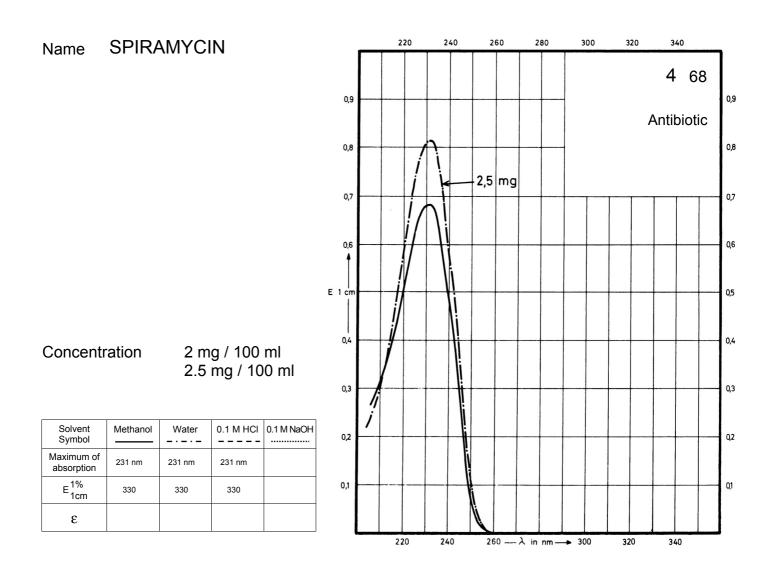


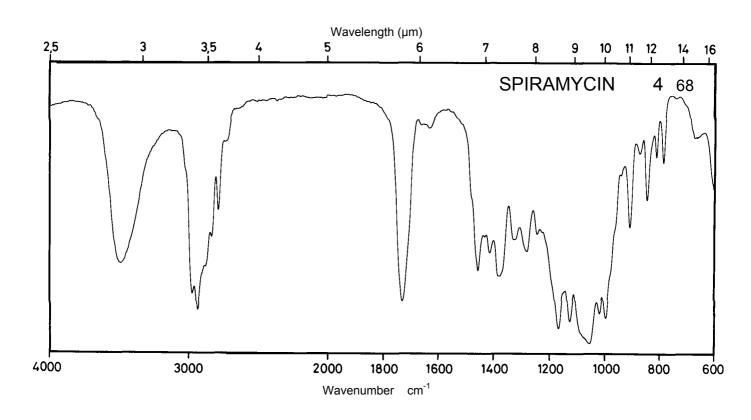

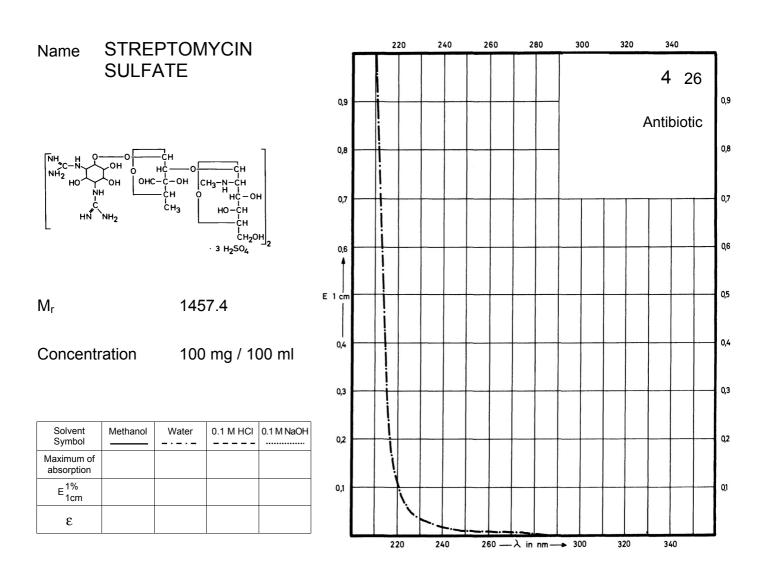

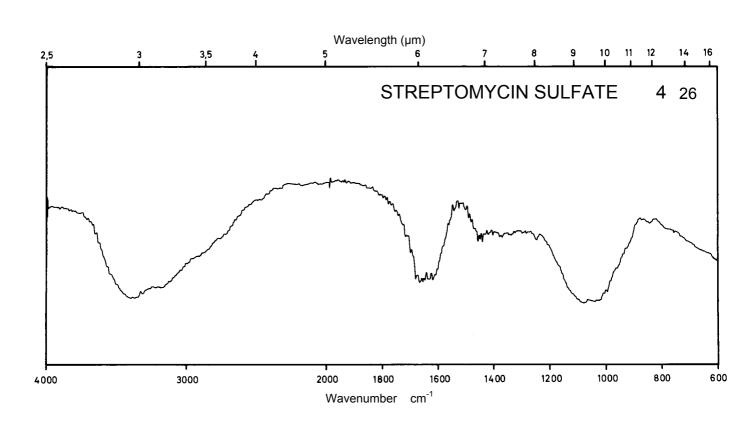


Name α -SANTONIN

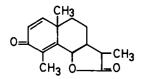


 M_{r}


246.3

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	238 nm		244 nm	247 nm
E ^{1%} 1cm	505		541	470
ε	12440		13320	11580

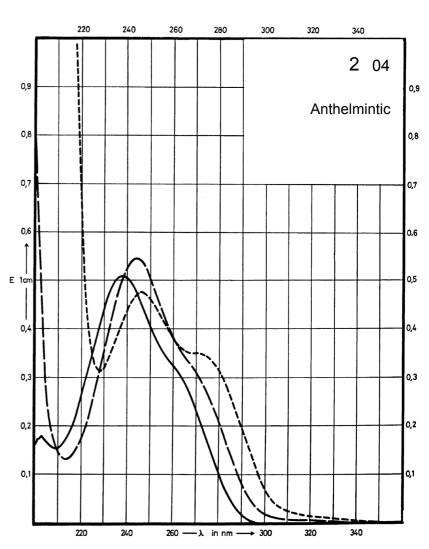

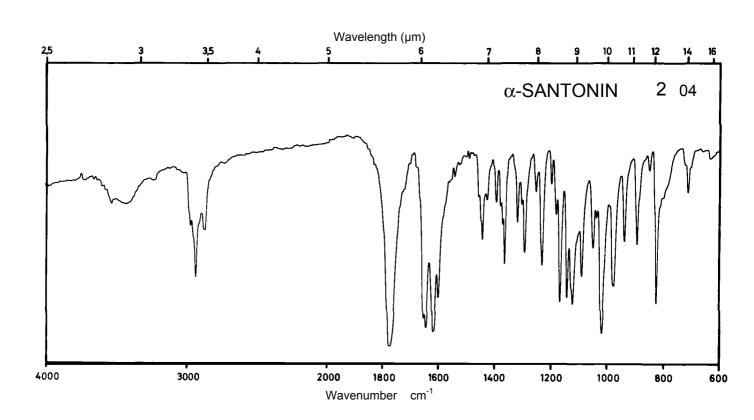



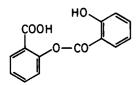



258.2

Concentration


1.6 mg / 100 ml

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	308 nm 235 nm		307 nm 237 nm	340 nm
E ^{1%} 1cm	190 587		177 659	254
3	4900 15160		4570 17000	6560

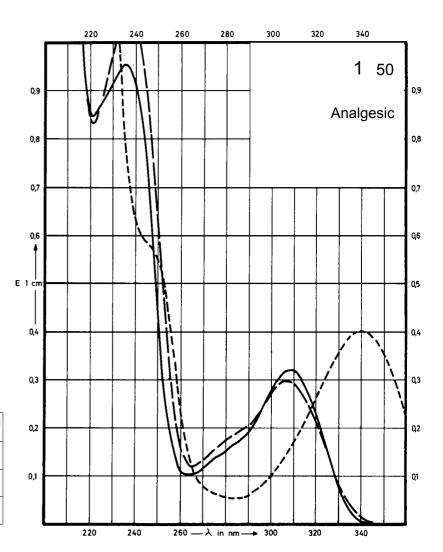

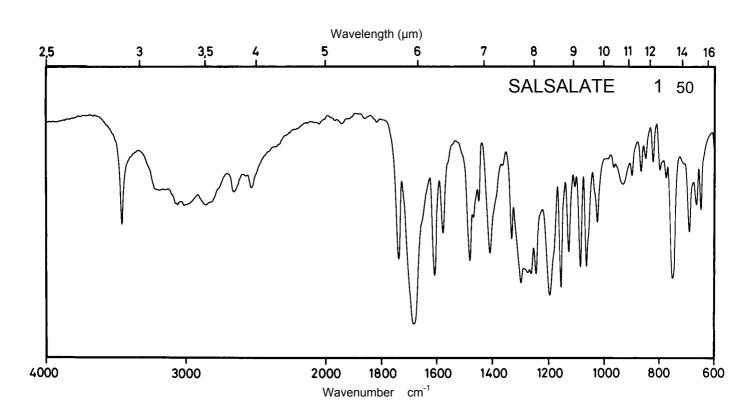





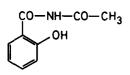


 M_r 179.2

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	307 nm 243 nm		308 nm 248 nm	Decom- position observed
E ^{1%} _{1cm}	229 710		211 658	
3	4100 12700		3800 11800	

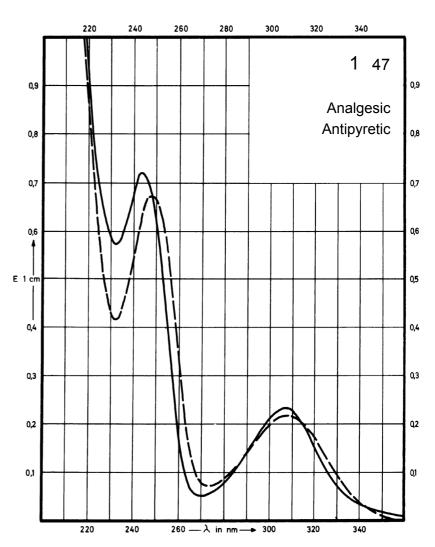

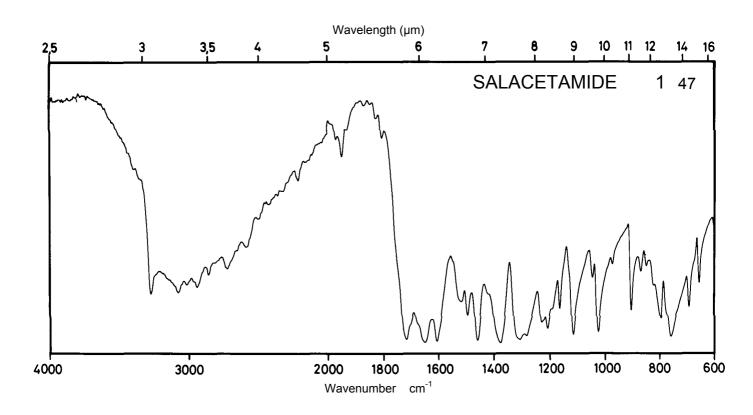





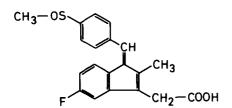

 M_r 356.4

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	328 nm 287 nm 226 nm		327 nm 284 nm 228 nm	327 nm 280 nm
E ^{1%} 1cm	375 432 547		366 419 559	364 442
ε	13370 15380 19370		13050 14920 19900	12980 15770

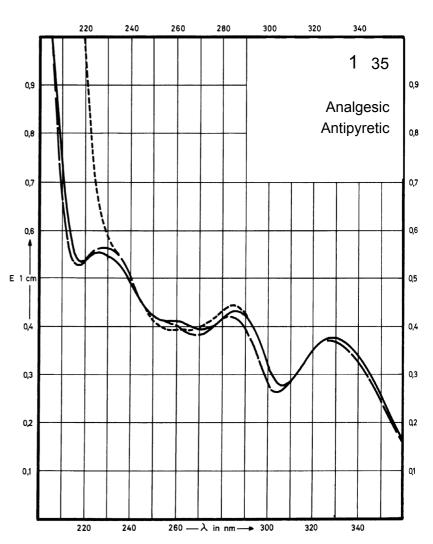

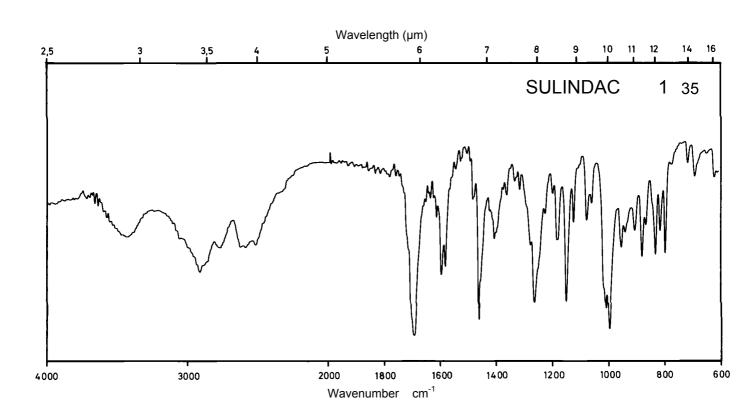







137.1

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	303 nm 235 nm		298 nm 234 nm	327 nm 240 nm
E ^{1%} 1cm	288 532		254 555	425 524
ε	3950 7290		3480 7610	5830 7180

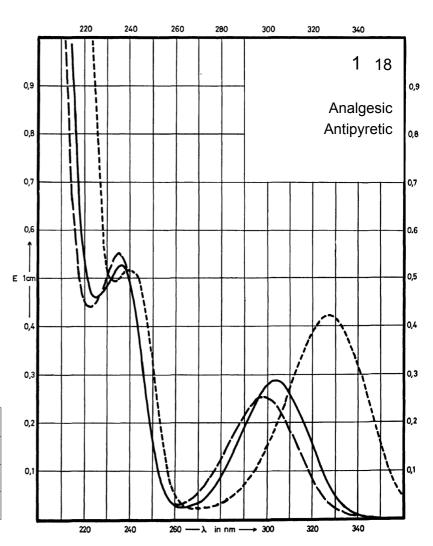

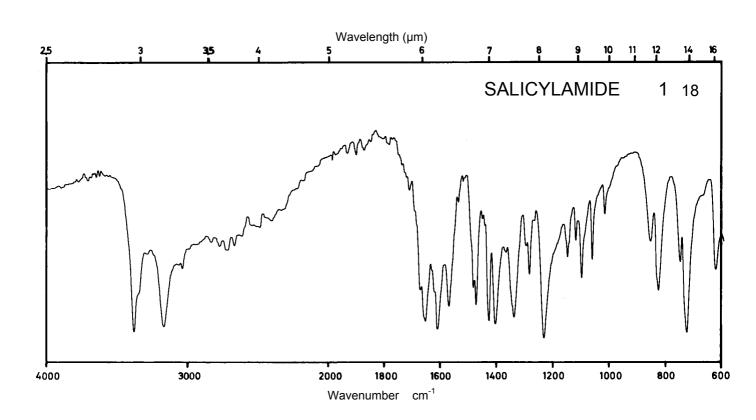



Name STRYCHNINE SULFATE

 M_{r}


857.0

Solvent Symbol	Methanol	Water	0.1 M HCI	0.1 M NaOH
Maximum of absorption	254 nm		254 nm	255 nm
E ^{1%} 1cm	300		298	291
ε	25700		25540	24940

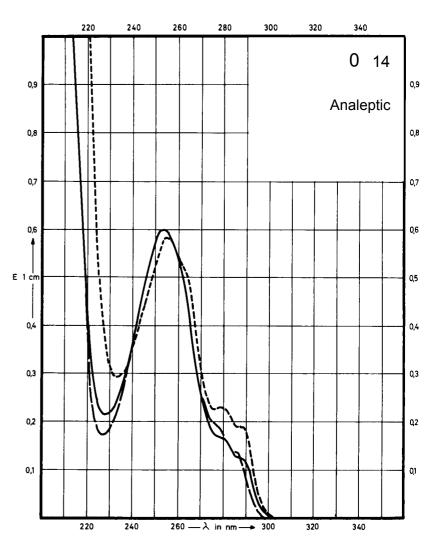


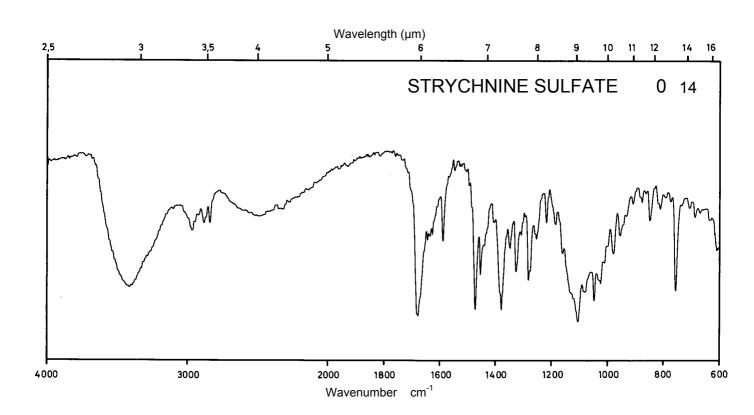


Name	Gelatine	Sample preparation Solid film, sodium chloride cell
		Reference
		Air
		Ordinate
		0 - 100 %

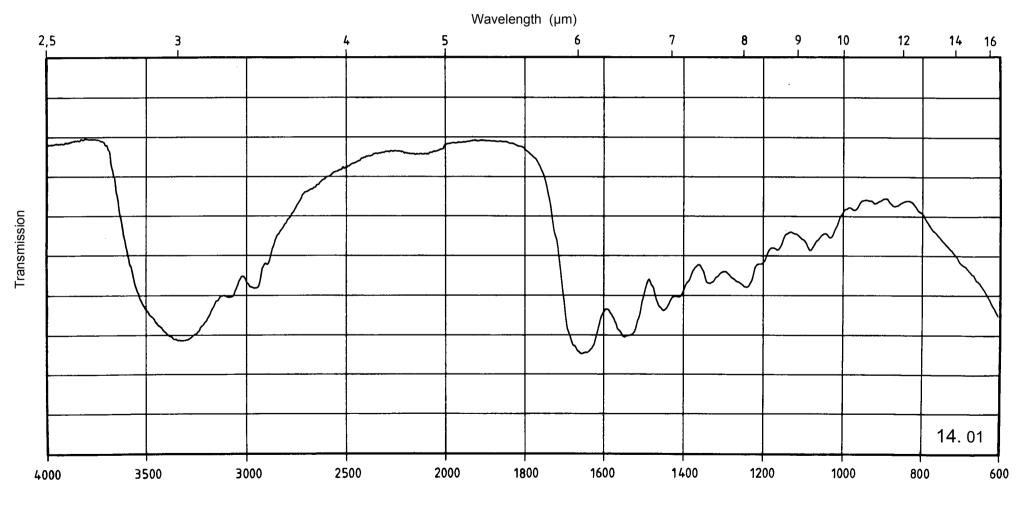


Name Propyl	ene glycol	Sample preparation Capillary film, sodium chloride cell
		Reference Air
		Ordinate 0 - 100 %

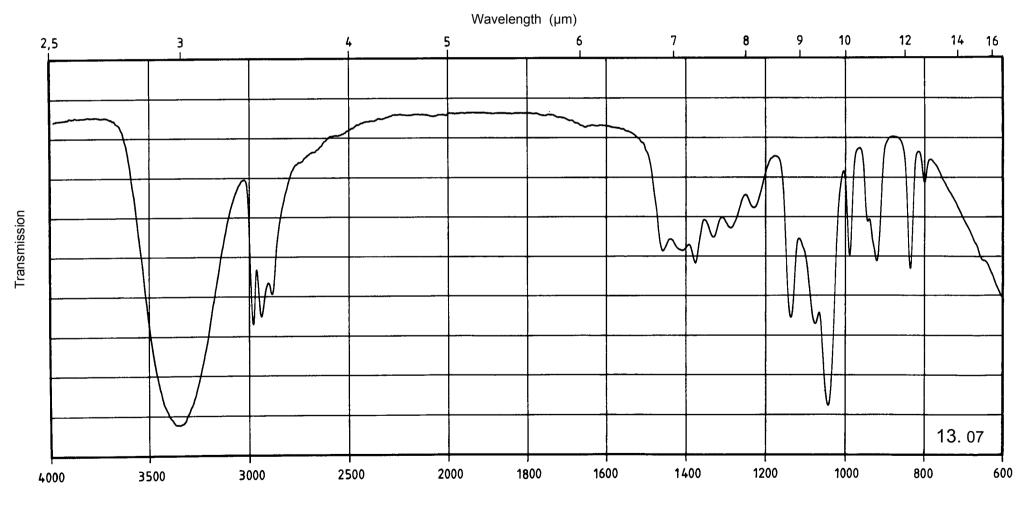


Name	Polyoxyethylene polypropylene glycol	Sample preparation Capillary film, sodium chloride cell
	Genapol® PF-10	Reference Air
		Ordinate 0 - 100 %

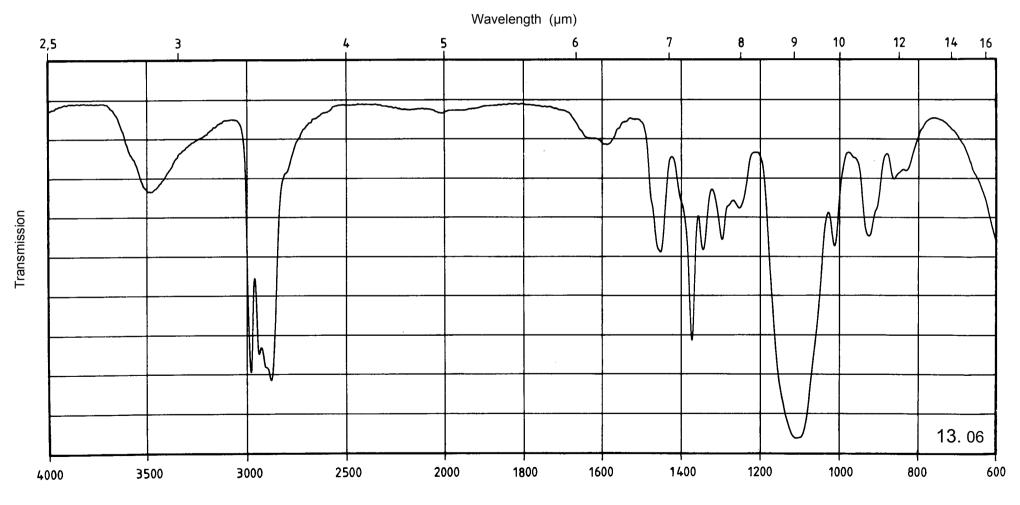


Name Polyethylene Glycol, PEG 4000	Sample preparation	
Nume	1 diyetilyiche diyedi, 1 20 1000	Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

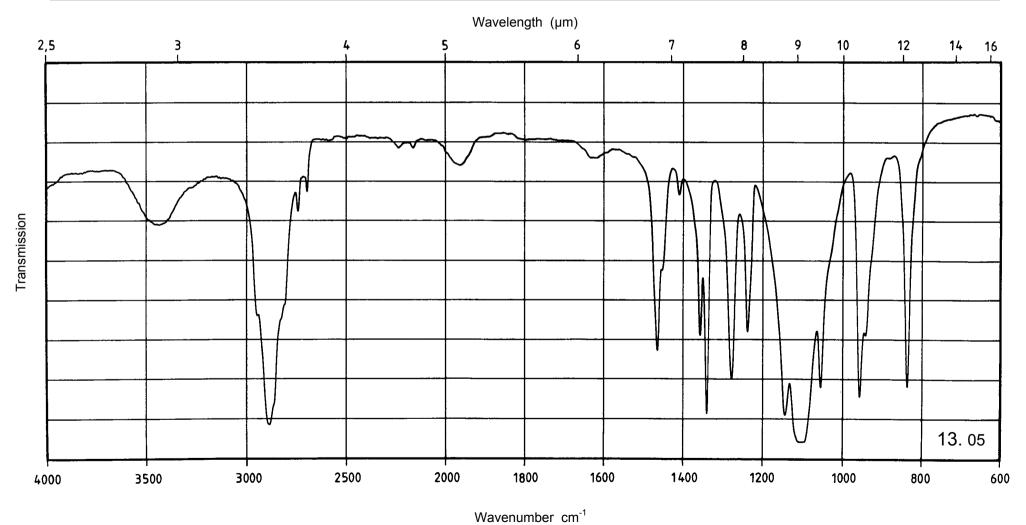


Name	Polyethylene glycol, PEG 400	Sample preparation Capillary film, sodium chloride cell
		Reference Air Ordinate 20 - 100 %

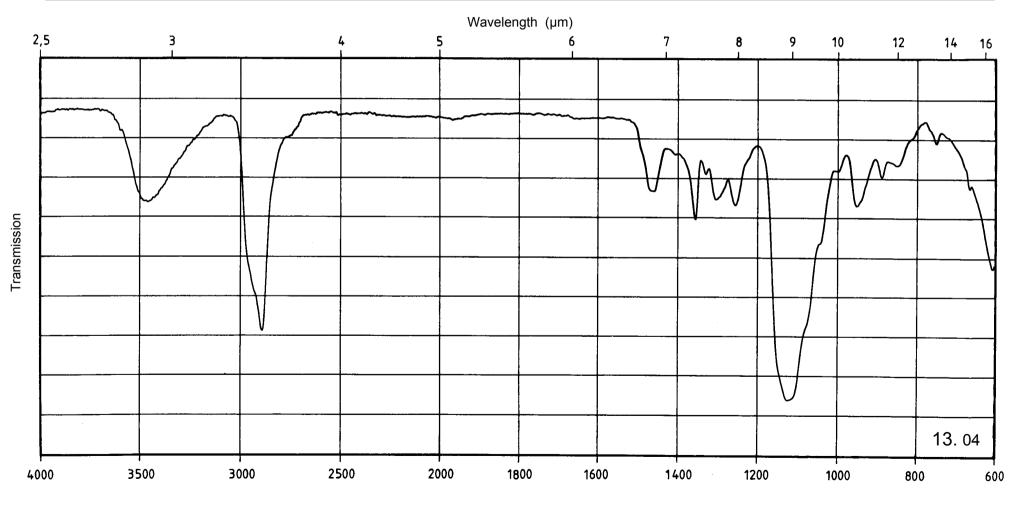



Name	Glycerol formal	Sample preparation Capillary film, sodium chloride cell
		Reference
		Air
		Ordinate
		0 - 100 %

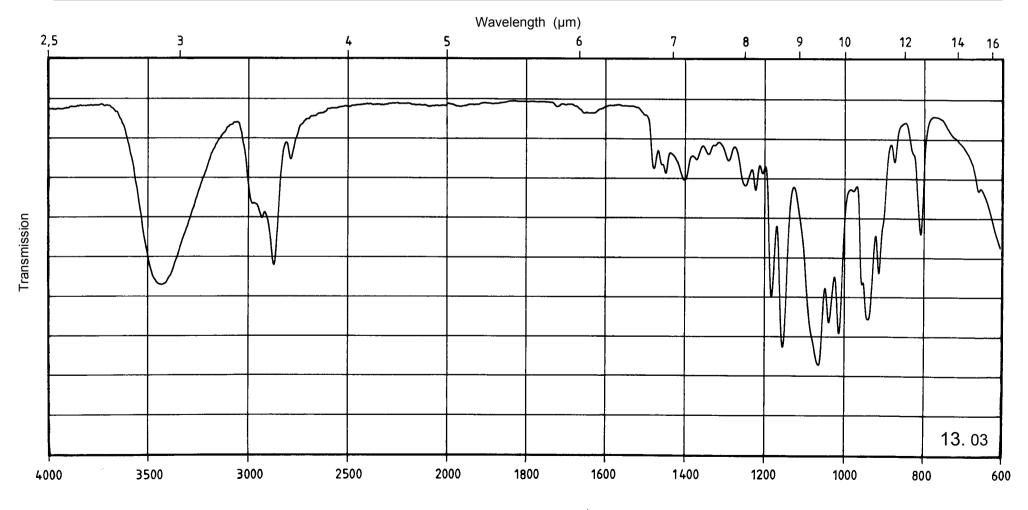


Name	Glycerol, 88 per cent	Sample preparation Capillary film, sodium chloride cell
		Reference
		Air
		Ordinate
		0 - 100 %

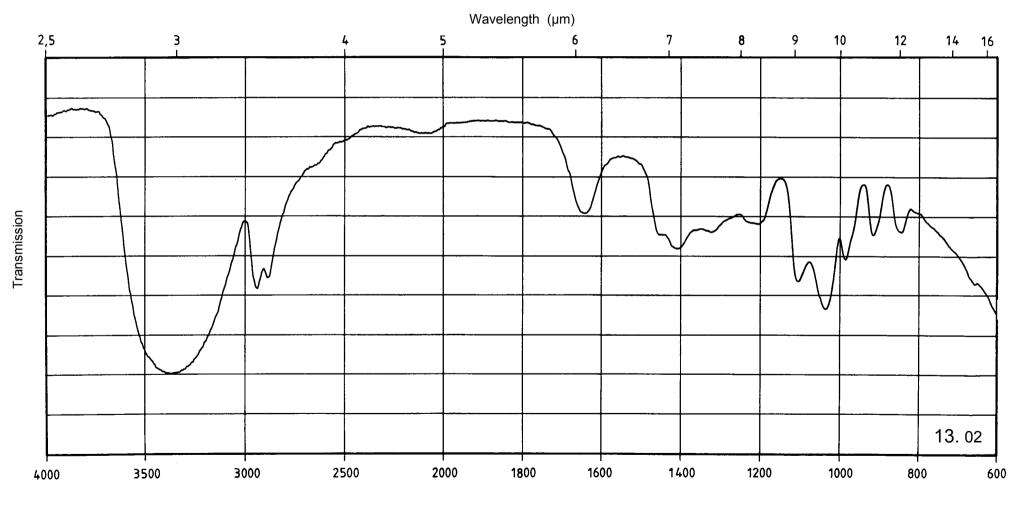

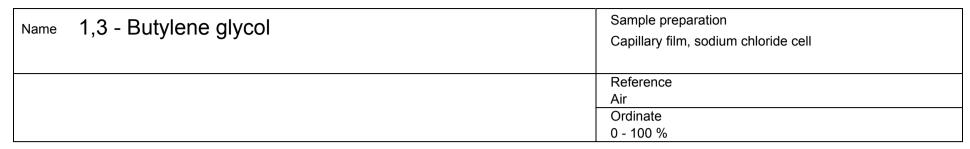


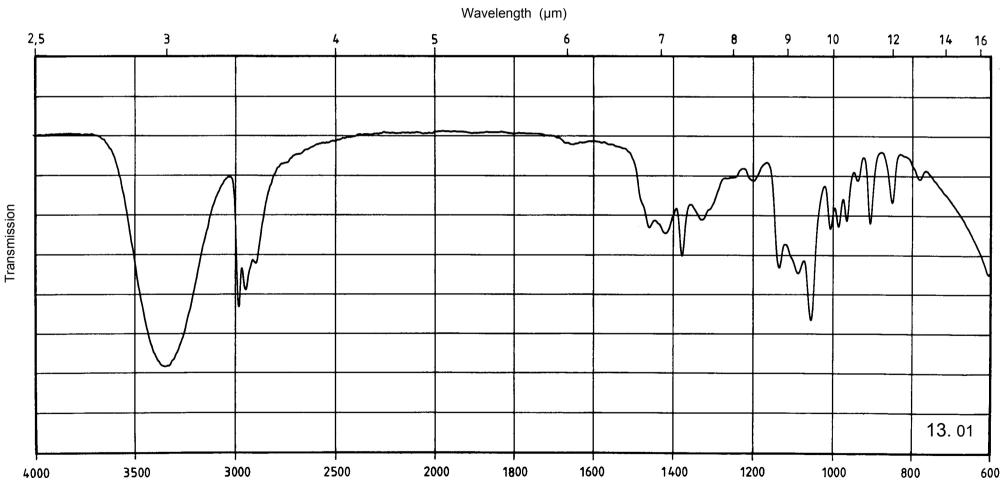


Name	Sorbitan, tristearate	Sample preparation
		Potassium bromide dispersion
		1.5 mg / 300 mg
	Span [®] 65	Reference
	Span 00	Air
		Ordinate
		0 - 100 %

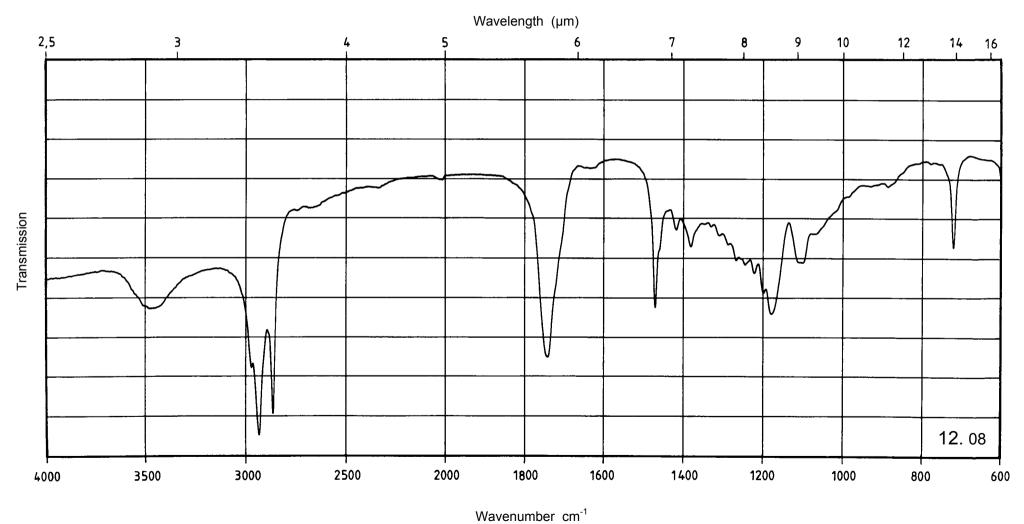


Name	Sorbitan, monostearate	Sample preparation Potassium bromide dispersion
		0.7 mg / 300 mg
	Span® 60	Reference Air
		Ordinate 0 - 100 %

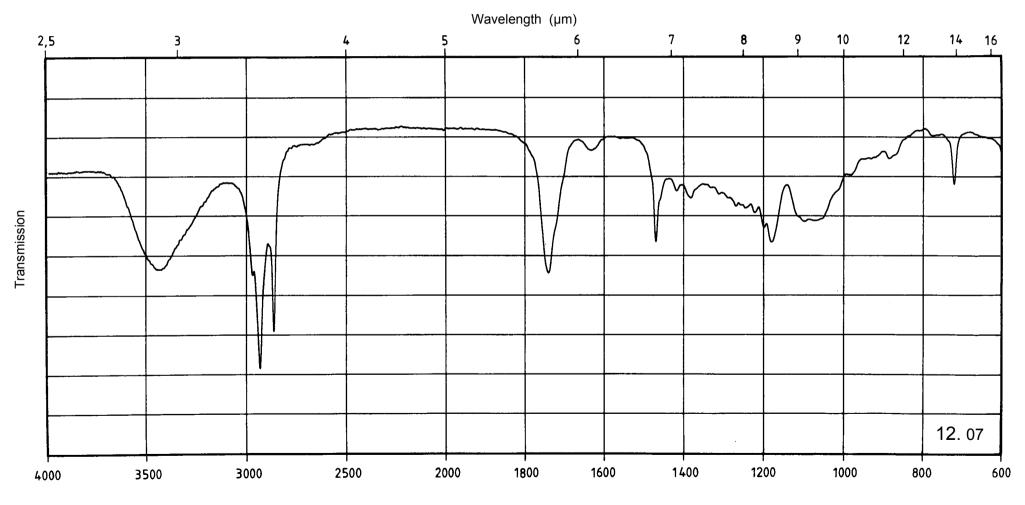


Name	Sorbitan, monooleate	Sample preparation Capillary film, sodium chloride cell
	Span [®] 80	Reference Air
		Ordinate 0 - 100 %

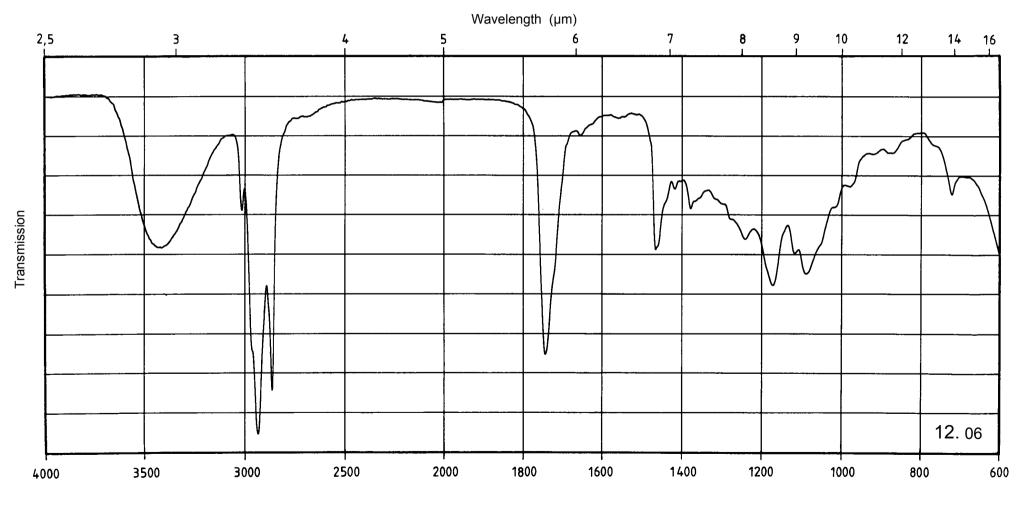


Name	Polyoxyethylene(20) sorbitan trioleate	Sample preparation Capillary film, sodium chloride cell
	Tween [®] 85, Polysorbat 85	Reference Air
		Ordinate 0 - 100 %

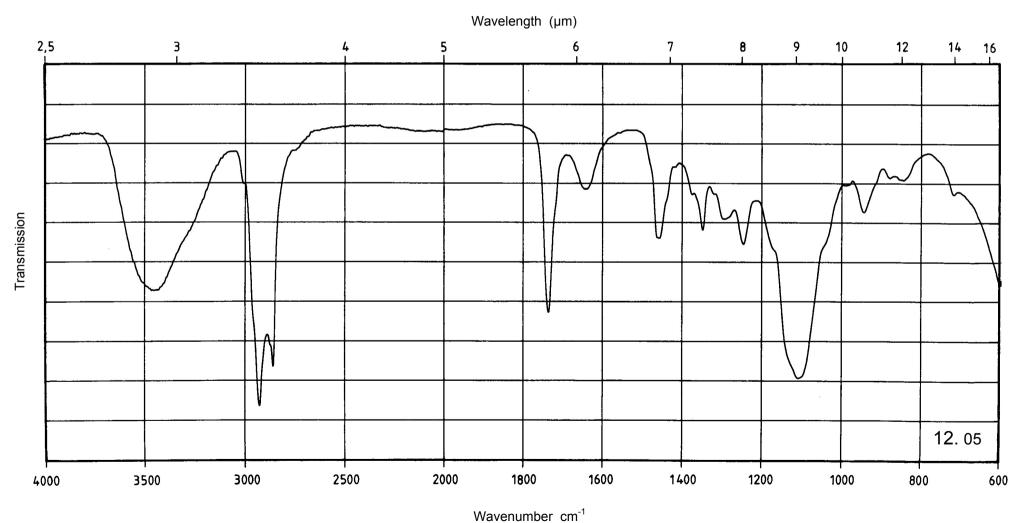


Name	Polyoxyethylene(20) sorbitan monostearate	Sample preparation Capillary film, sodium chloride cell
	Tween [®] 60, Polysorbat 60	Reference Air
		Ordinate 0 - 100 %

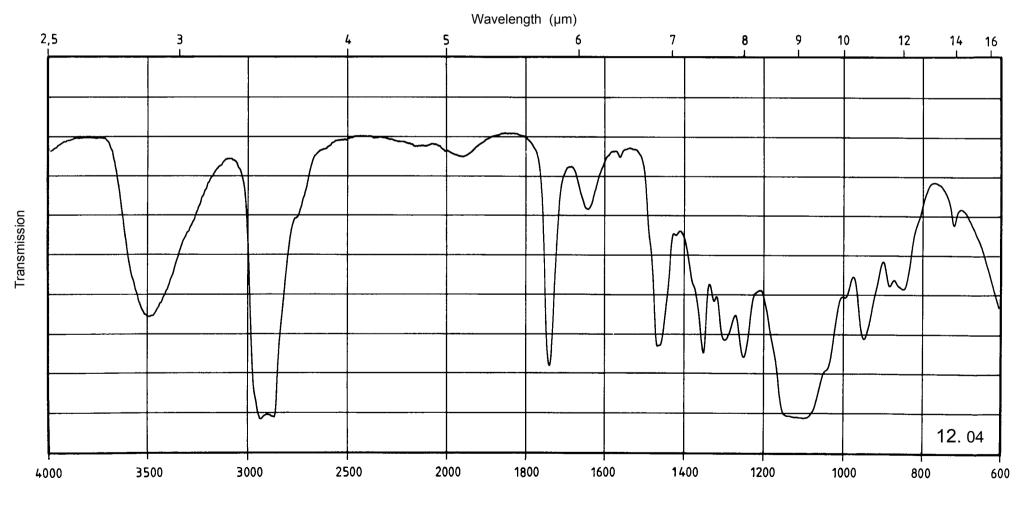


Name	Polyoxyethylene(20) sorbitan monooleate	Sample preparation Capillary film, sodium chloride cell
	Tween [®] 80, Polysorbat 80	Reference Air Ordinate
		0 - 100 %

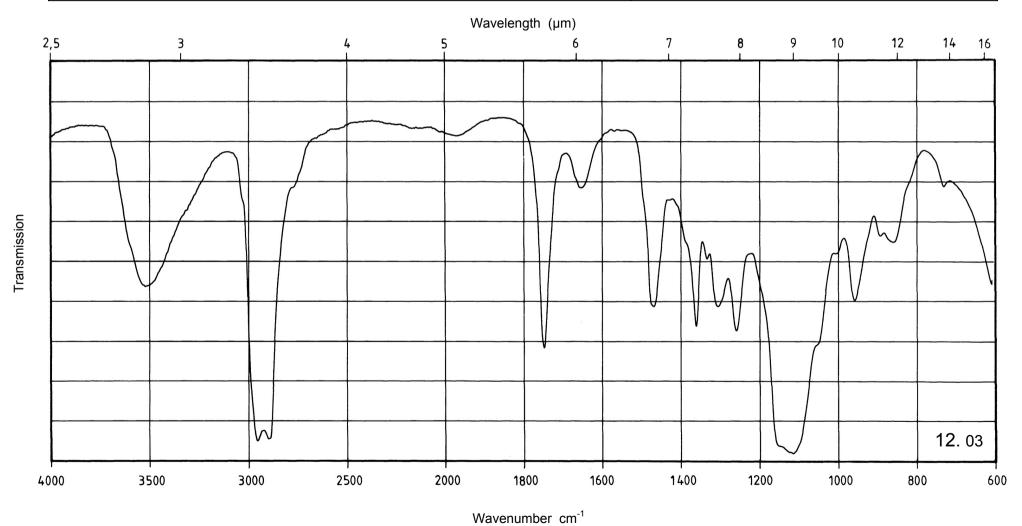


Name	Polyoxyethylene(20) sorbitan monolaurate	Sample preparation Capillary film, sodium chloride cell
	Tween [®] 20, Polysorbat 20	Reference Air
		Ordinate 0 - 100 %

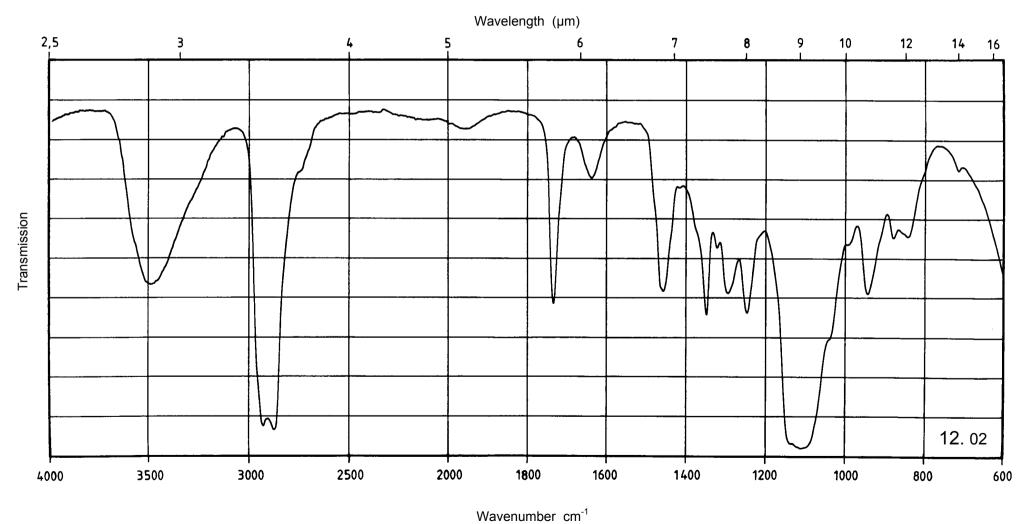


Name	Glycerol sorbitan stearate	Sample preparation Capillary film, sodium chloride cell
	Arlacel® 481	Reference Air
		Ordinate 0 - 100 %

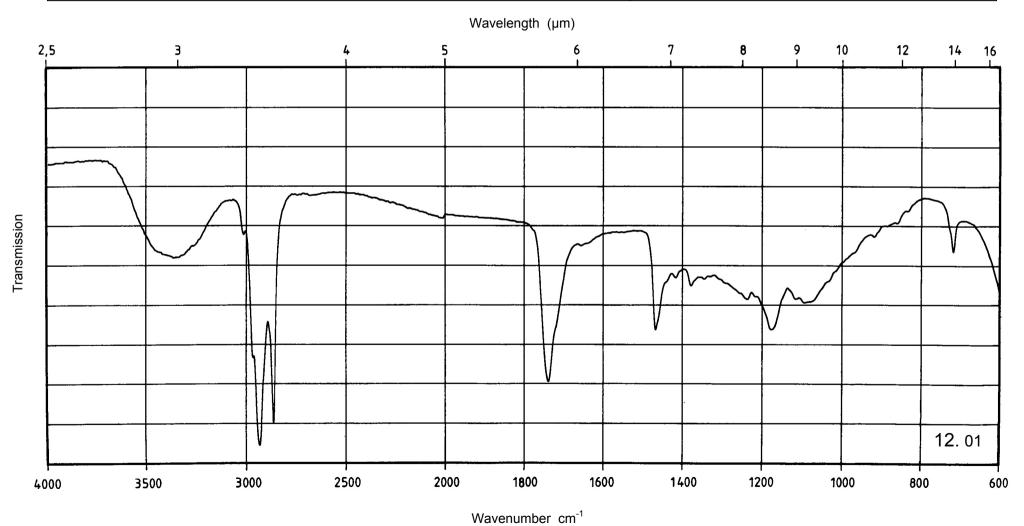


Name	1-Vinyl-2-pyrrolidone-Polymer, cross-linked Crosspovidone	Sample preparation Potassium bromide dispersion
	Orosopovidorio	1.5 mg / 300 mg
	Kallidan Cl. Dalynlaedana VI	Reference
	Kollidon CL, Polyplasdone XL	Air
		Ordinate
		0 - 100 %

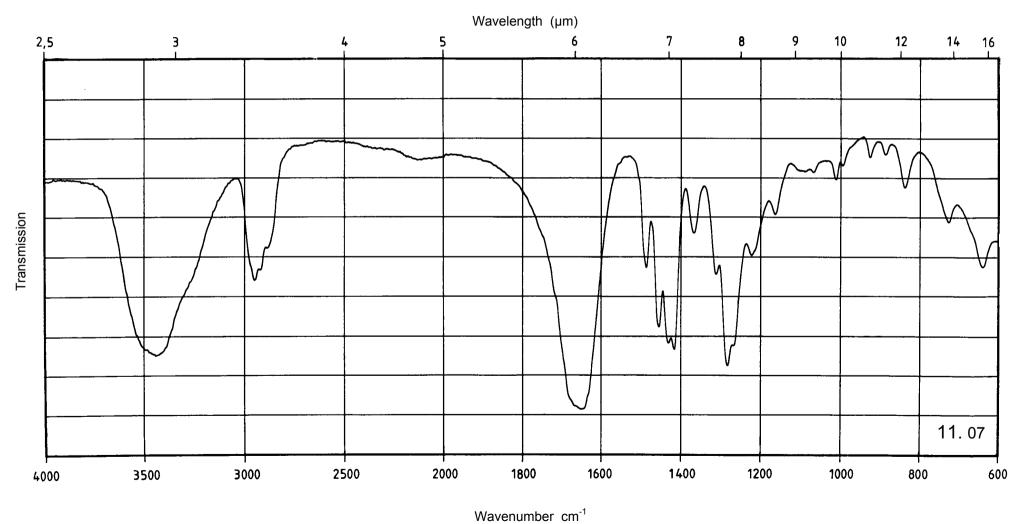


Name 1-Vinyl-2-pyrroli	1-Vinyl-2-pyrrolidone-Polymer	Sample preparation
	1 Villy 2 pyriolidenc i diyinci	Potassium bromide dispersion
		4 mg / 300 mg
	Kollidon ® 30. Povidone	Reference
	Kollidon 30. Povidone	Air
		Ordinate
		0 - 100 %

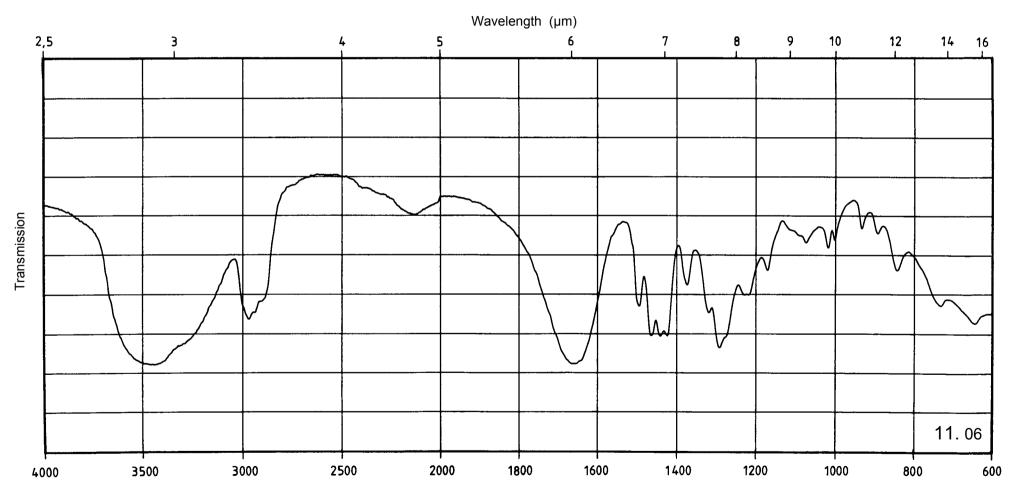


Name	Polyvinyl acetate	Sample preparation Solid film, sodium chloride disk
		Reference
		Air
		Ordinate
		0 - 100 %

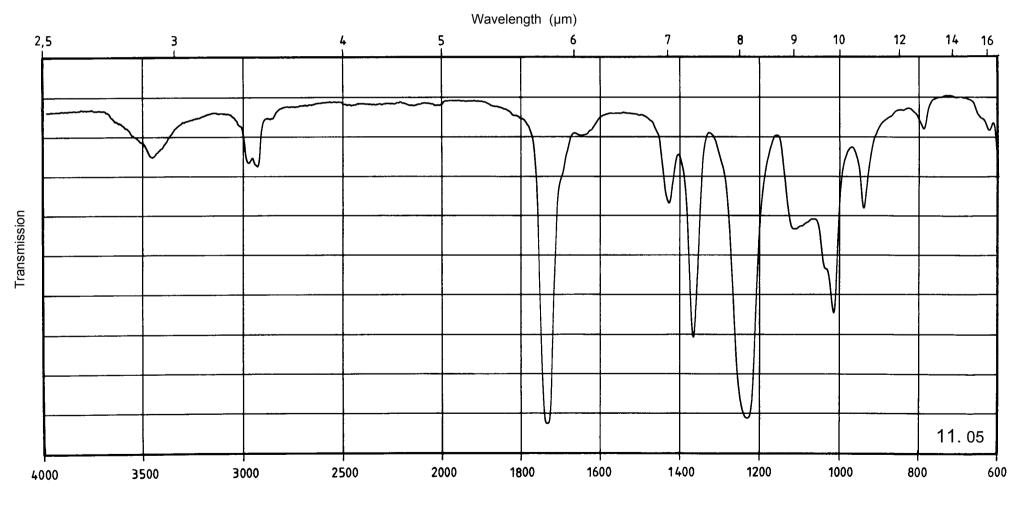


Name	Polyisobutene	Sample preparation Solid film, sodium chloride disk
	Oppanol B100	Reference Air
		Ordinate
		0 - 100 %

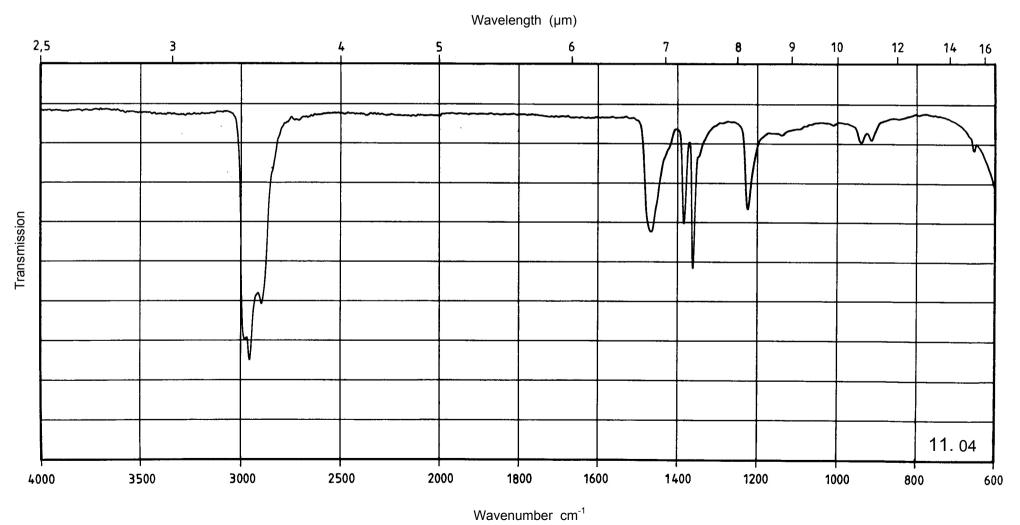


Name	Poly-3-hydroxybutyric acid	Sample preparation
	1 ory o riyaroxybatyric dola	Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

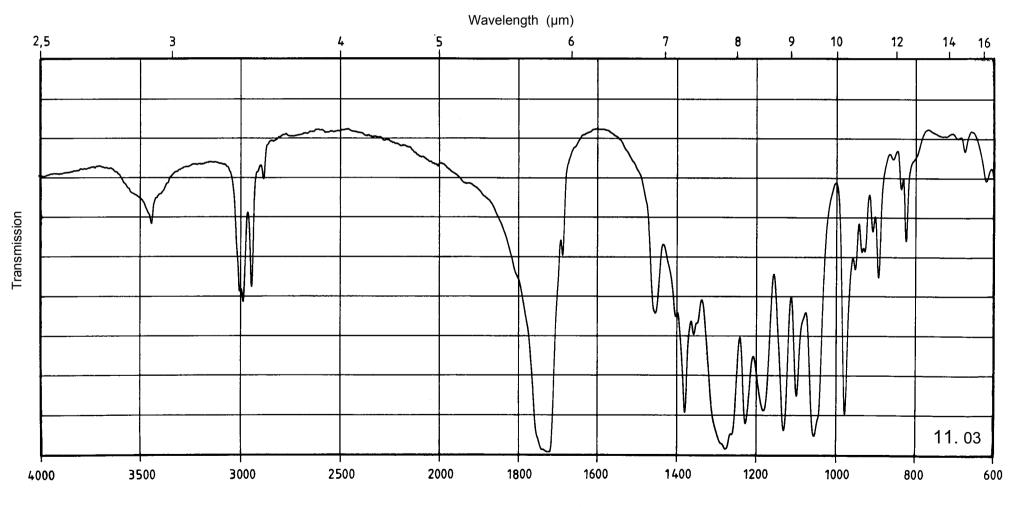


Name	Carboxypolymethylene	Sample preparation
	Carboxyporymetrylene	Potassium bromide dispersion
		1.5 mg / 300 mg
	Carbopol 941 [®]	Reference
	Carbopol 94 i	Air
		Ordinate
		0 - 100 %

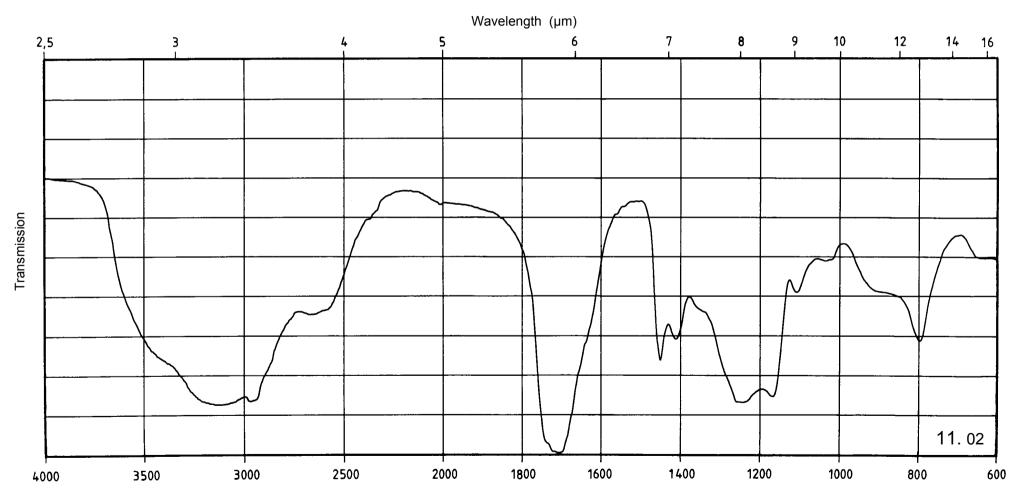


Name	Name Methacrylic acid, methacrylic acid ester, copolymer	Sample preparation
Ivanie		Potassium bromide dispersion
		1.5 mg / 300 mg
	Eudragit ® L100	Reference
	Eudragit L100	Air
		Ordinate
		0 - 100 %

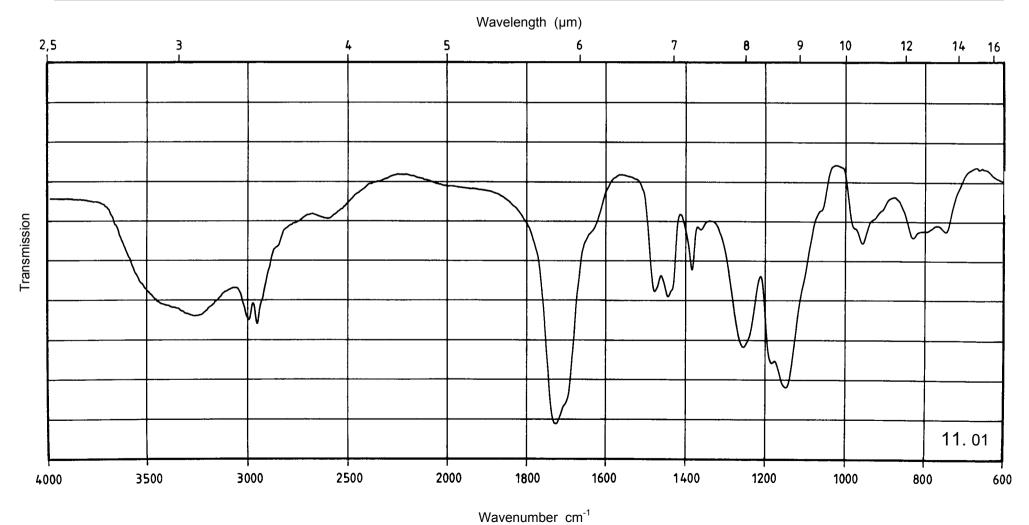


Name	Phenylmethylpolysiloxane	Sample preparation Capillary film, sodium chloride cell
		Reference
		Air Ordinate 0 - 100 %

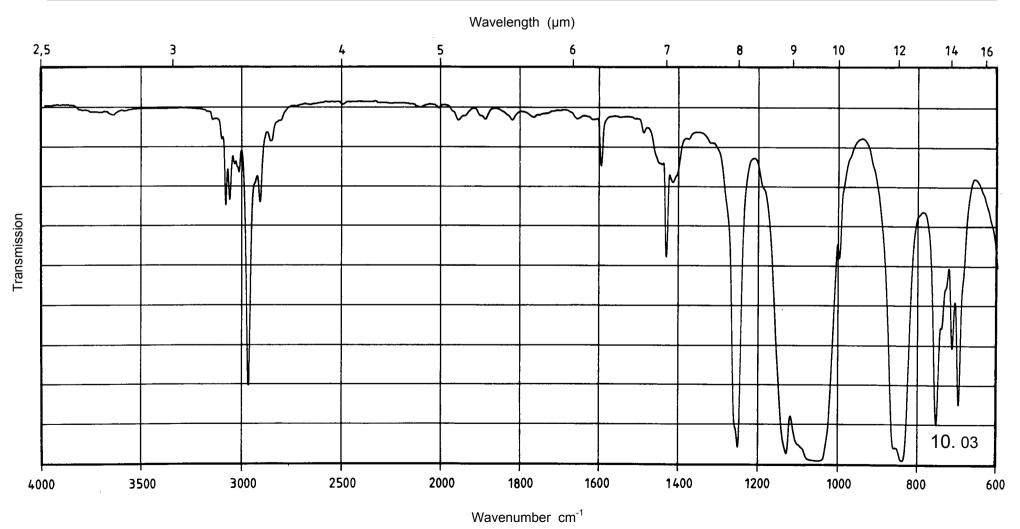


Name	Polydimethylcyclosiloxane	Sample preparation Capillary film, sodium chloride cell
		Reference Air
		Ordinate 0 - 100 %

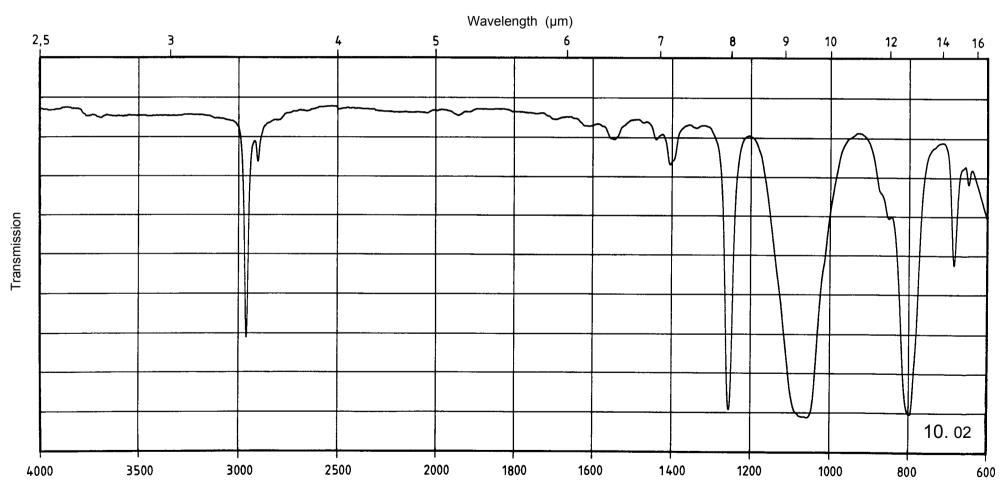


Name	Dimethylpolysiloxane	Sample preparation Capillary film, sodium chloride cell
		Reference Air
		Ordinate 0 - 100 %

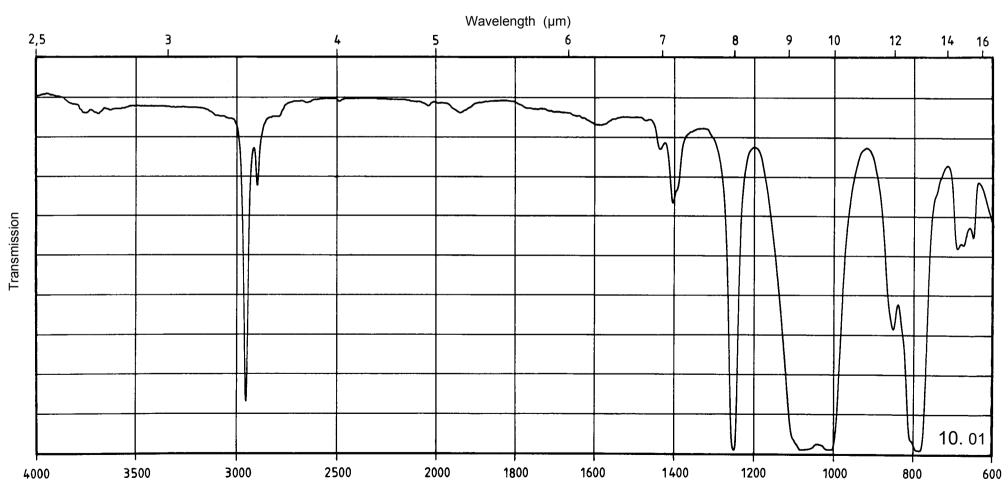


Name Talc	Sample preparation
Haine Paro	Potassium bromide dispersion
	1.5 mg / 300 mg
	Reference
	Air
	Ordinate
	0 - 100 %

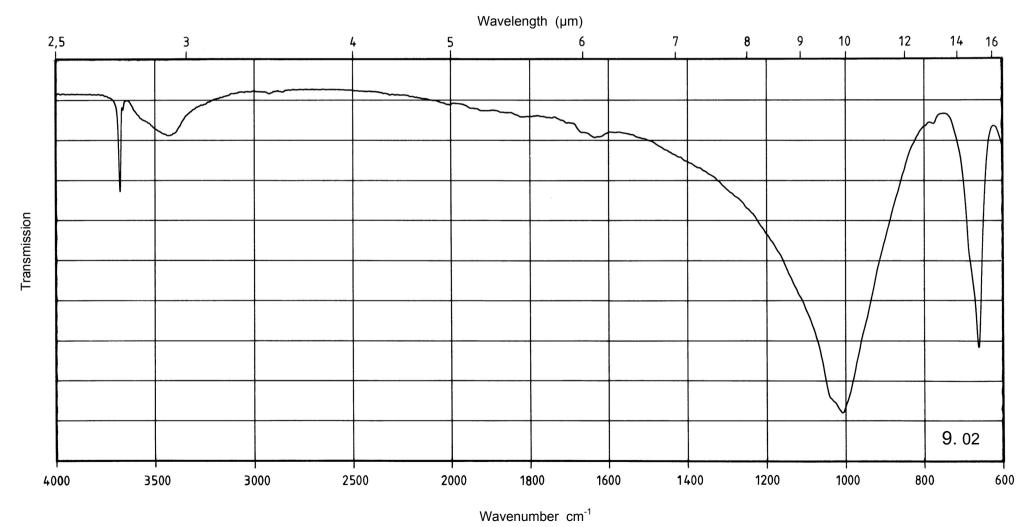


Name	Silicagel	Sample preparation
Ivallie	Ollicagei	Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

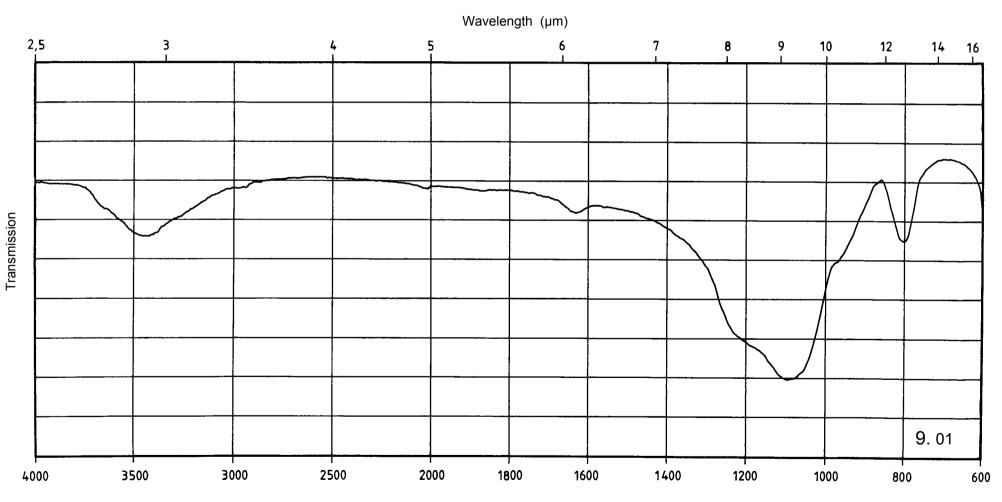


Name	Guar gum	Sample preparation Potassium bromide dispersion 3 mg / 300 mg
T Carrie	Caar gam	
		3 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

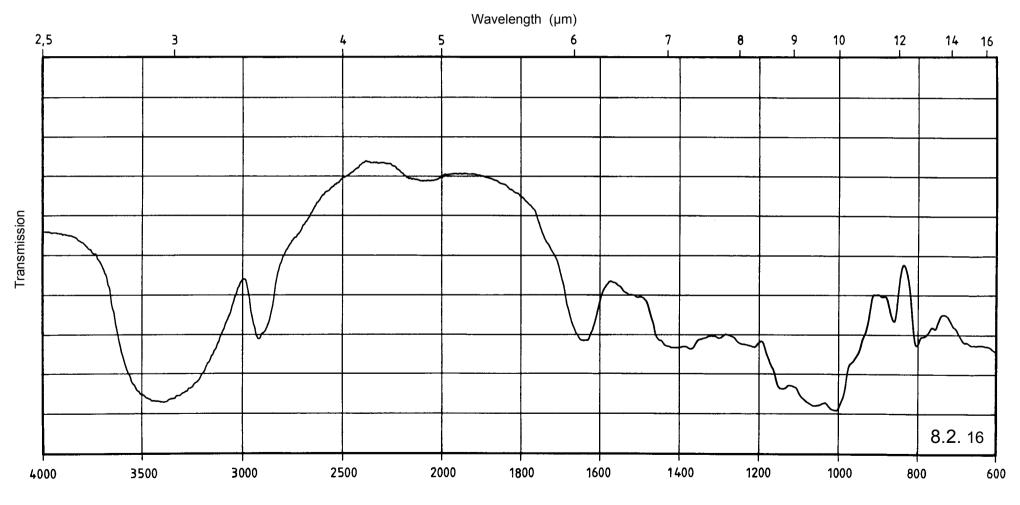


Name Hydroxypropylcellulose	Sample preparation Solid film
	Reference Air Ordinate 0 - 100 %

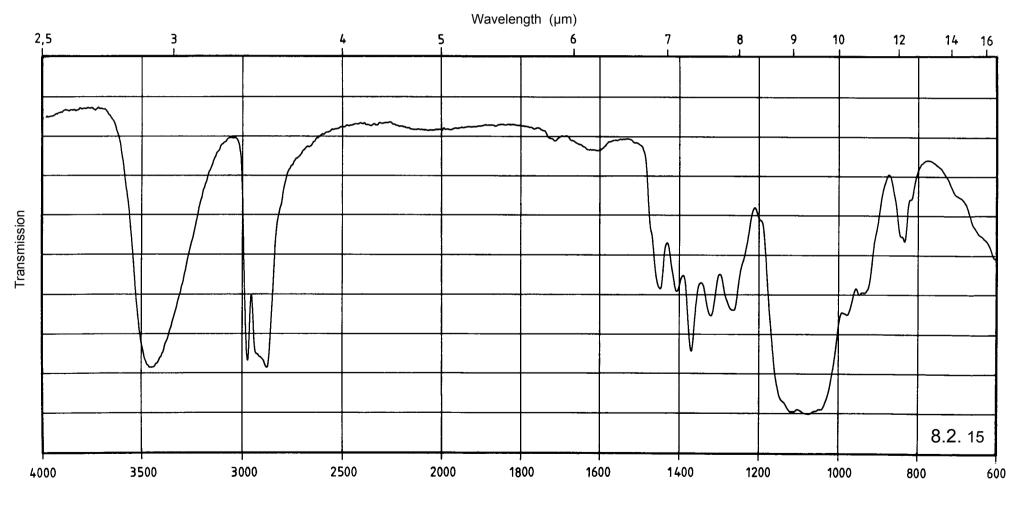


Name Hydroxyethylcellulose	Sample preparation Solid film
	Reference Air Ordinate 0 - 100 %

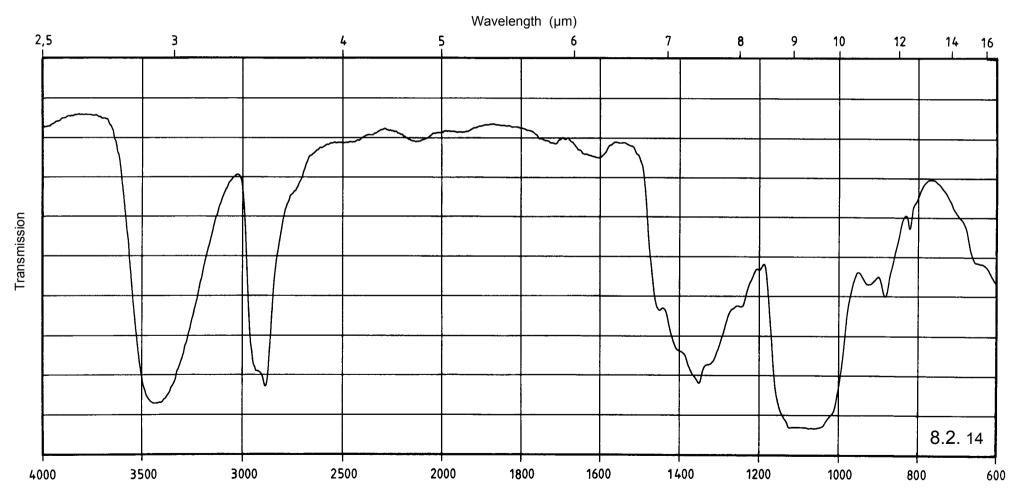


Name	Tragacanth	Sample preparation
		Potassium bromide dispersion
		3 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

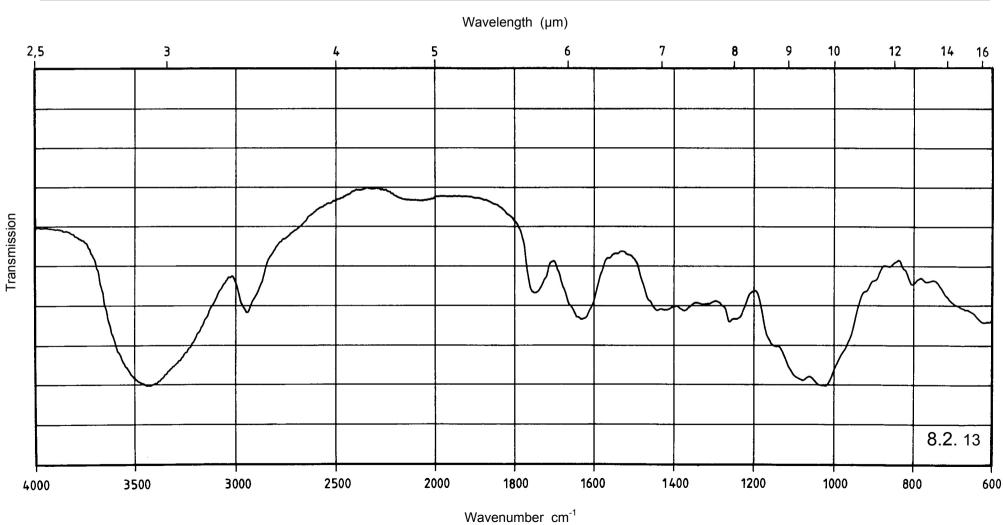


Name So	Sodium starch glycolate	Sample preparation
		Potassium bromide dispersion
		4 mg / 300 mg
	Primojel [®]	Reference
	Phinojei	Air
		Ordinate
		10 - 70 %

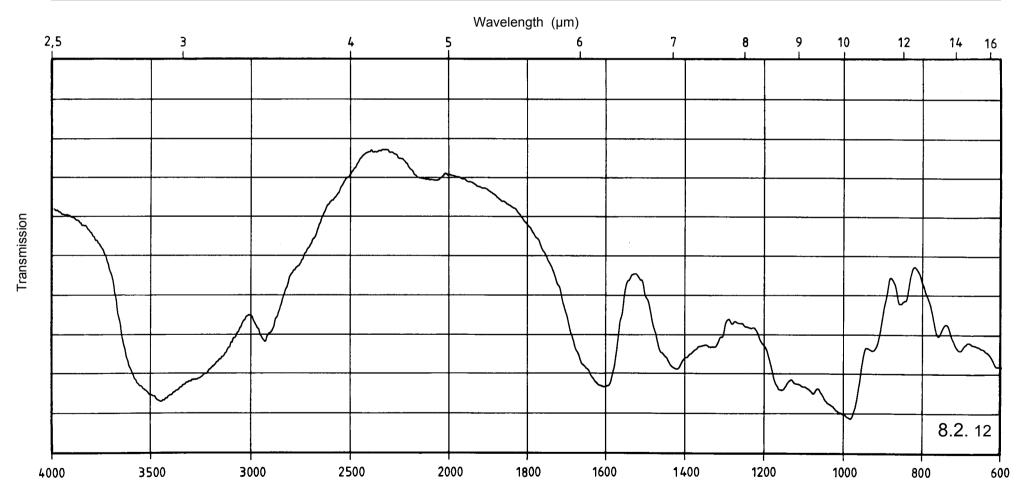


Name	Name Sodium carboxymethylcellulose	Sample preparation
Name		Potassium bromide dispersion
		4 mg / 300 mg
	Tylose® C 1000 P	Reference
	Tylose C 1000 P	Air
		Ordinate
		20 - 80 %

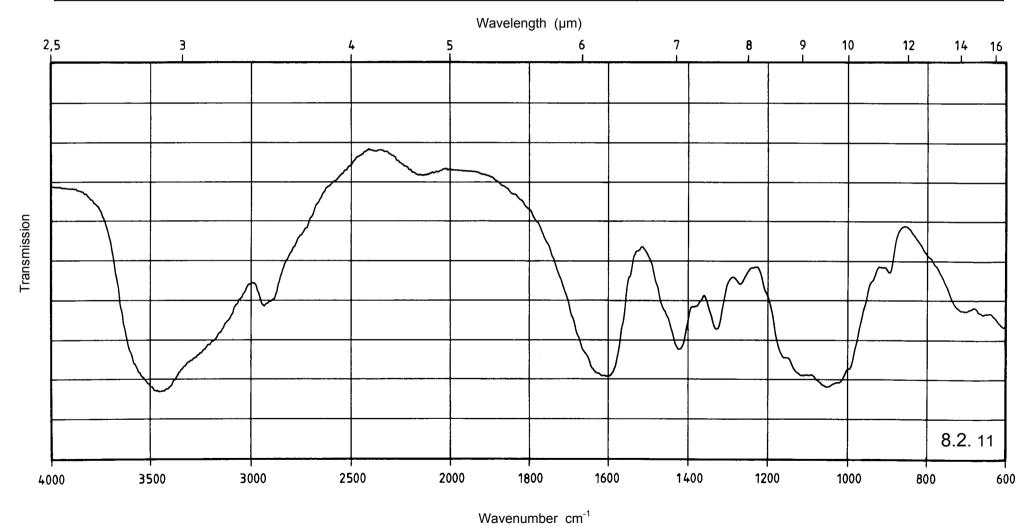


Name	Sodium amylopektin glycolate	Sample preparation
		Potassium bromide dispersion
		4 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

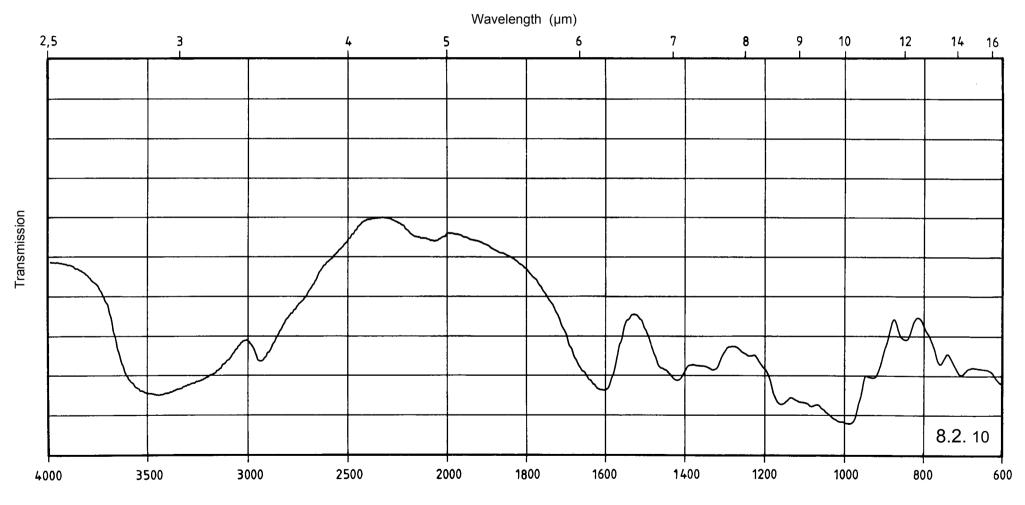


Name	Methyl hydroxyethylcellulose	Sample preparation Solid film, sodium chloride cell
	Tylose [®] MH 50	Reference Air
		Ordinate 0 - 100 %

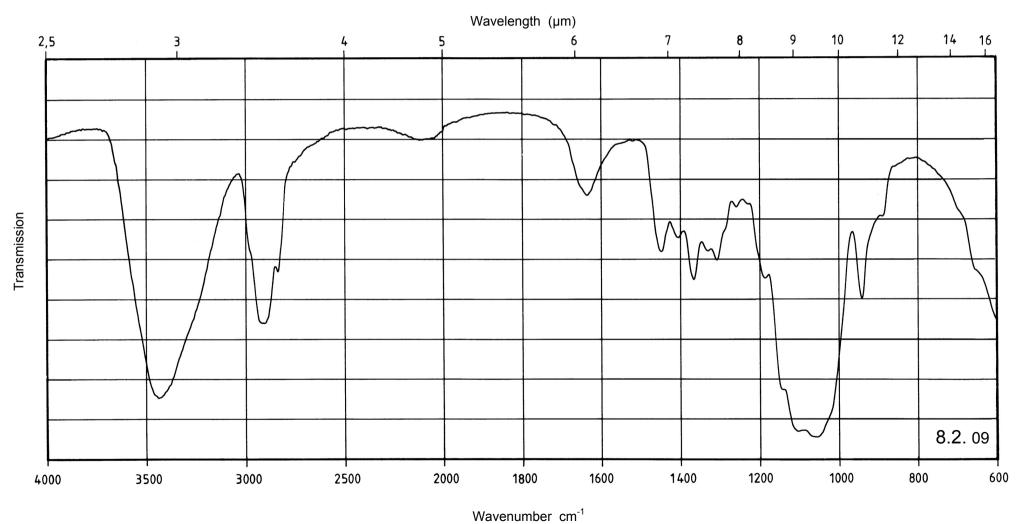


Name	Corn starch	Sample preparation Potassium bromide dispersion
		2 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

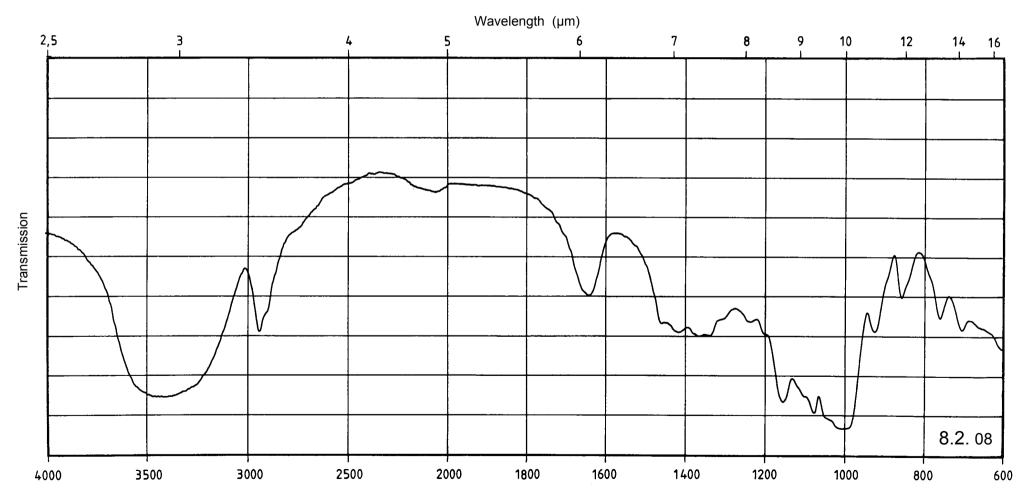


Name	Hydroxypropyl methyl cellulose phthalate	Sample preparation Solid film, potassium bromide disk
		Reference Air Ordinate 0 - 100 %

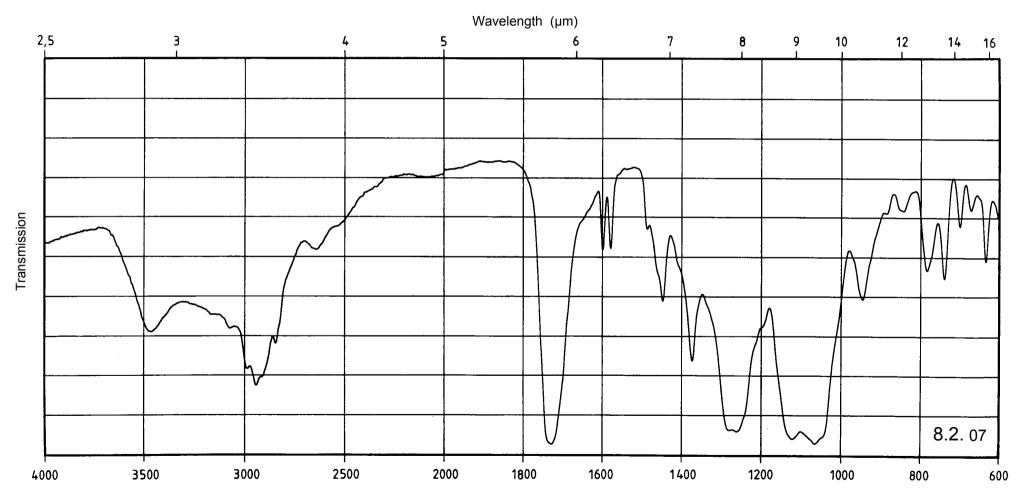


Name	Hydroxypropylmethylcellulose	Sample preparation
	Try di exty propy in round record	Potassium bromide dispersion
		4 mg / 300 mg
	Methocel [®] E5 Premium, Pharmacoat [®] 606	Reference
	Methodel E3 Premium, Pharmadoat 606	Air
		Ordinate
		0 - 100 %

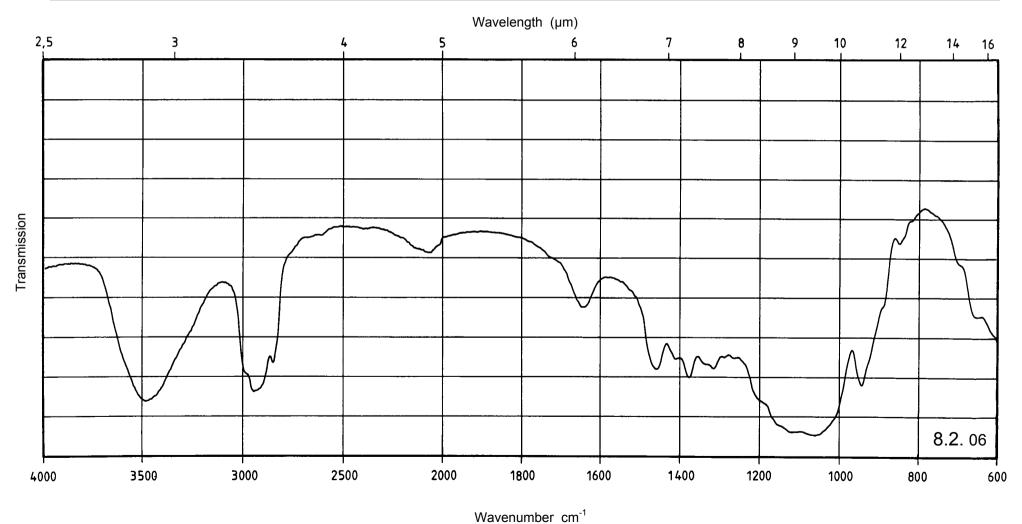


Name Acacia	Sample preparation
Name / todola	Potassium bromide dispersion
	3 mg / 300 mg
	Reference
	Air
	Ordinate
	0 - 100 %

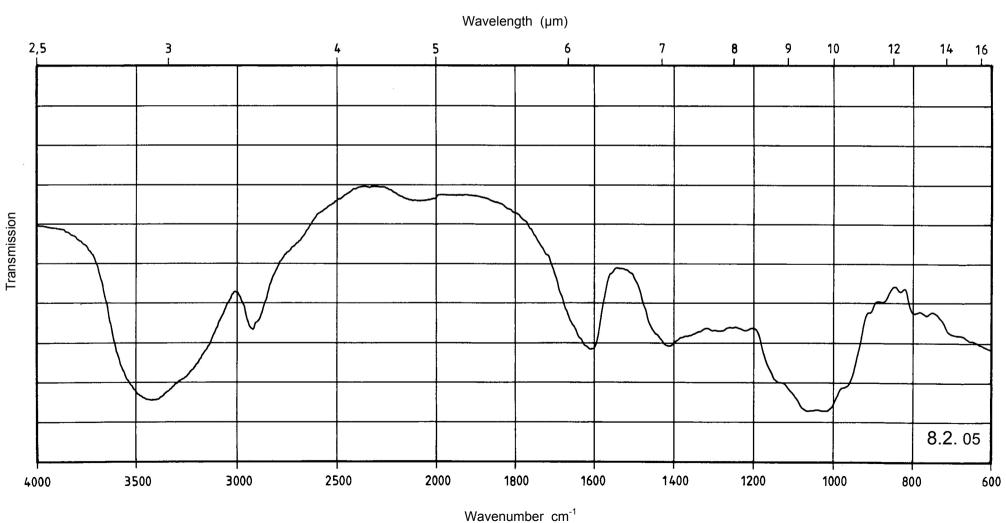


Name	Ethylcellulose	Sample preparation Solid film, potassium bromide disk
		Reference Air
		Ordinate 0 - 100 %

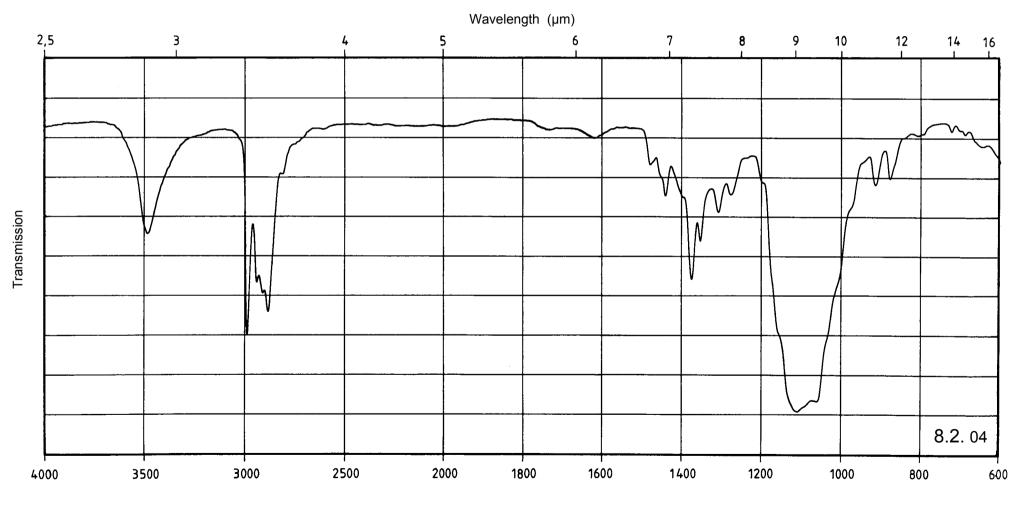


Name	β-Cyclodextrin	Sample preparation
		Potassium bromide dispersion
		4 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

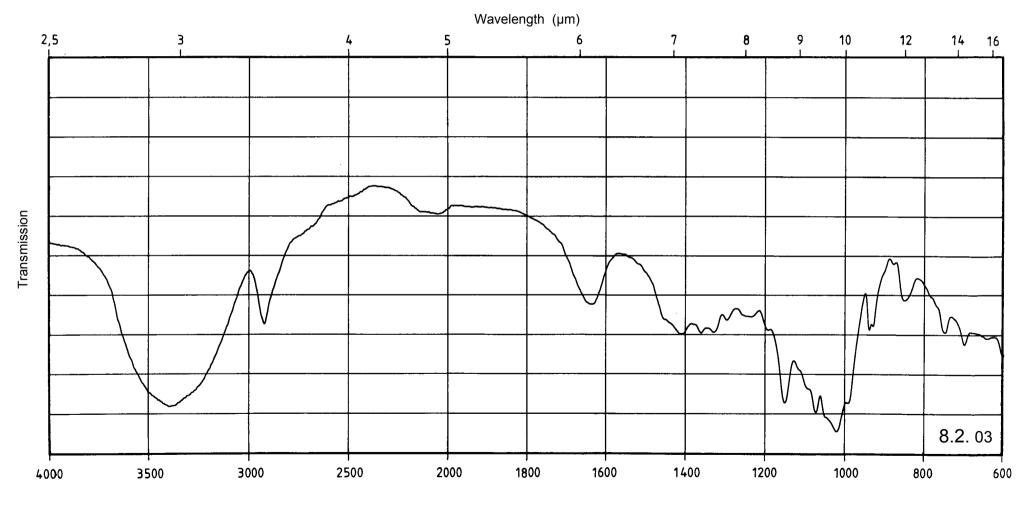


Name	Cellulose acetate phthalate	Sample preparation Solid film, potassium bromide disk
		Reference Air
		Ordinate 0 - 100 %

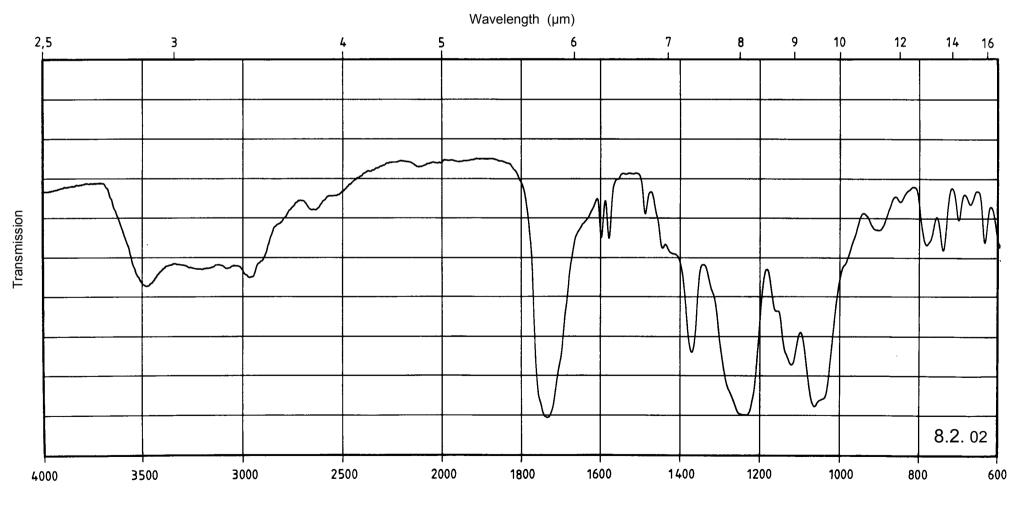


Name	e Cellulose	Sample preparation
Name		Potassium bromide dispersion
		3 mg / 300 mg
	Avicel® PH 100, Elcema® P100	Reference
	Avicei Ph 100, Elcella P100	Air
		Ordinate
		0 - 100 %

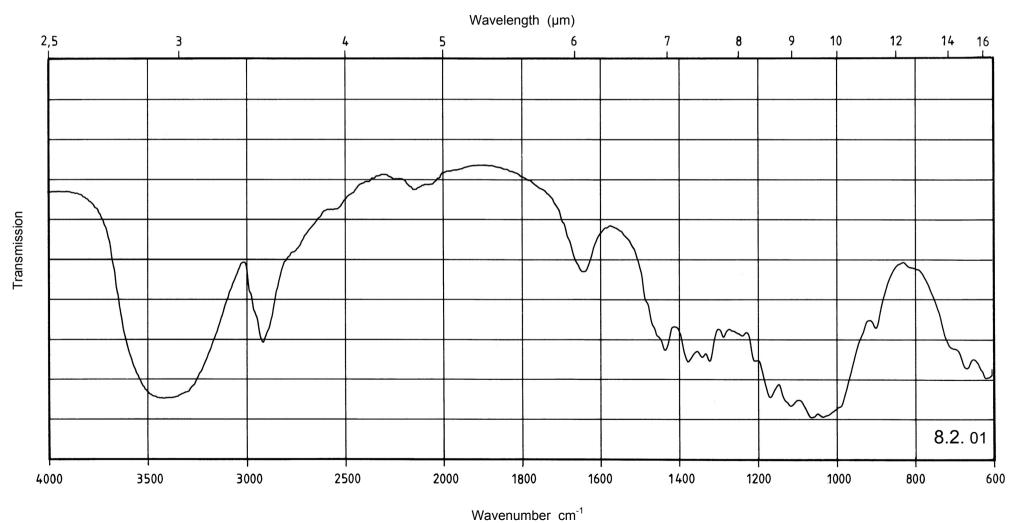


Name	e Sorbitol	Sample preparation
Ivanie		Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

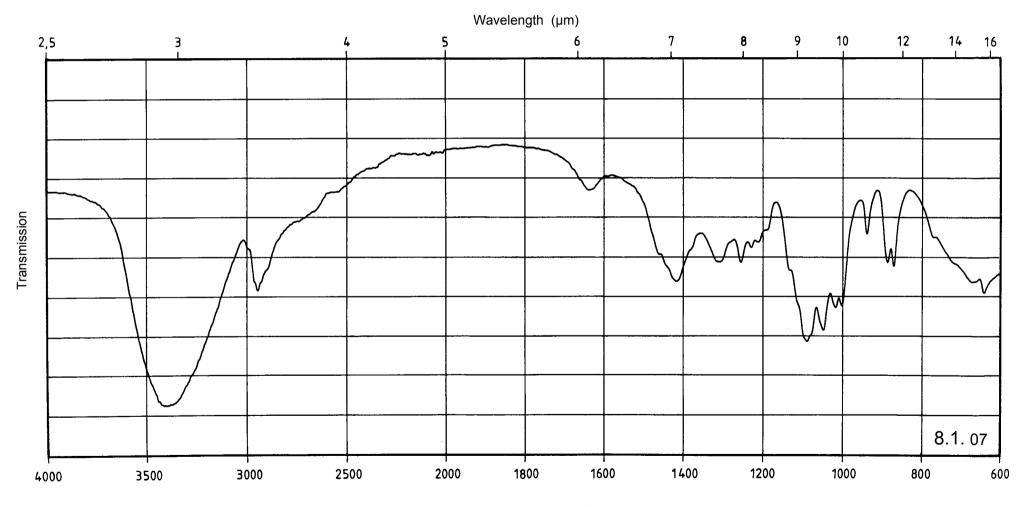


Name	Sucrose	Sample preparation
		Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

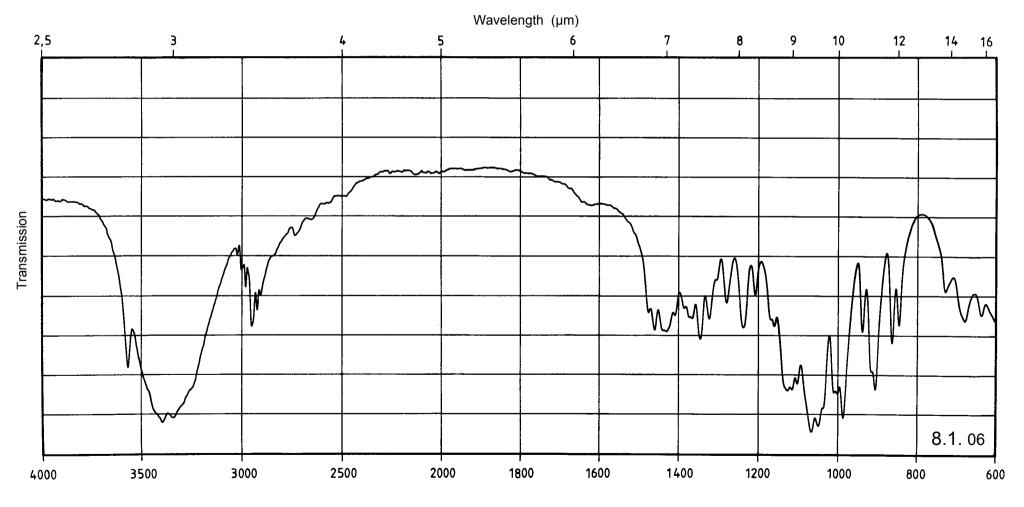


Name	e Mannitol	Sample preparation
Ivaille		Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

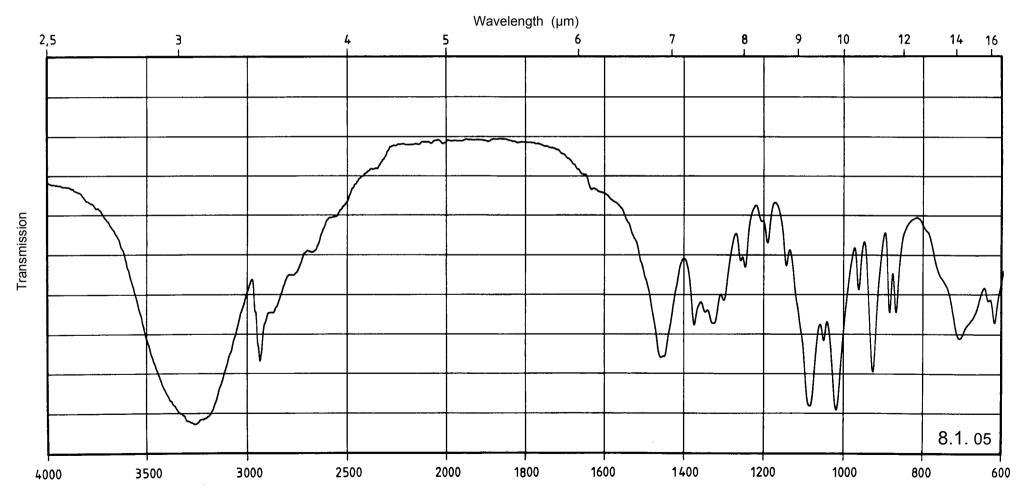


Name	Lactose	Sample preparation Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

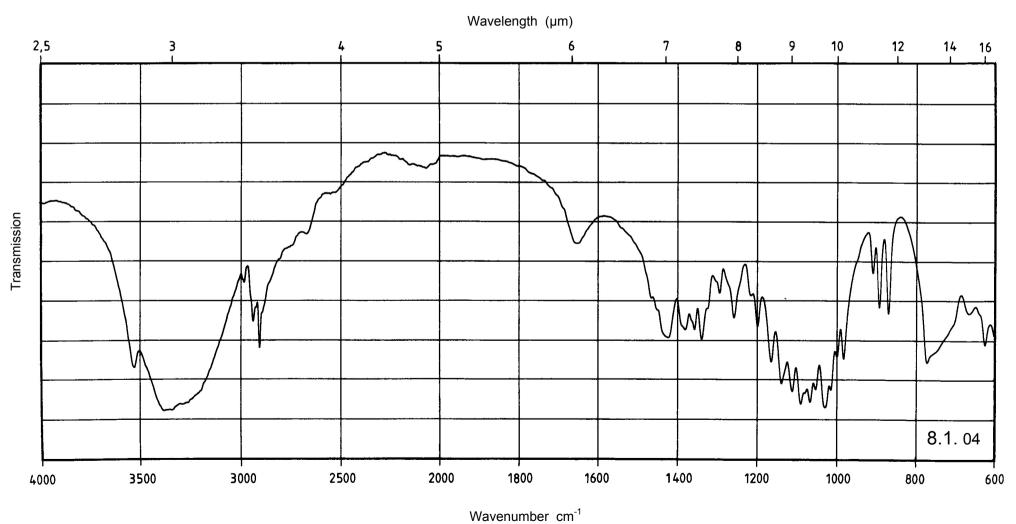


Name	Glucose, anhydrous	Sample preparation
Ivaille		Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

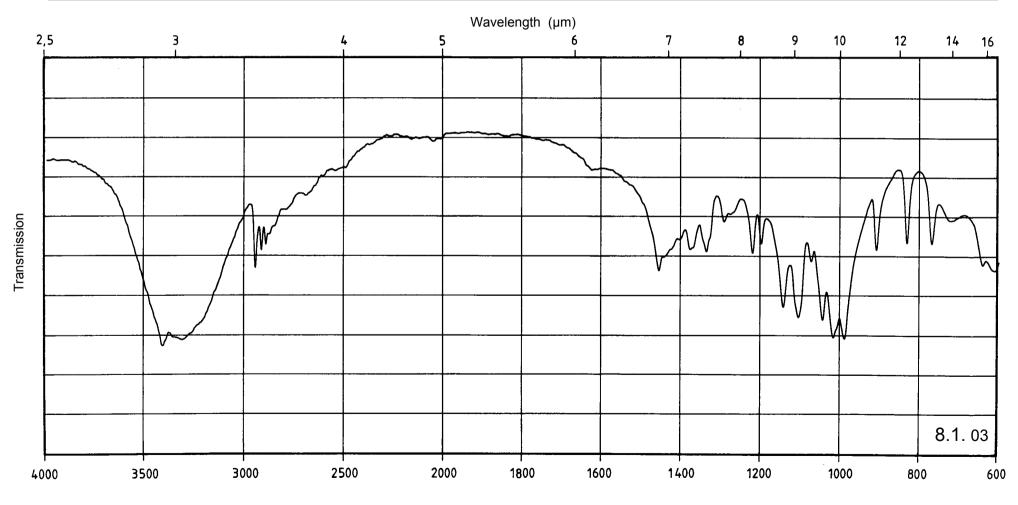


Name	Glucose, monohydrate	Sample preparation
		Potassium bromide dispersion
		2 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

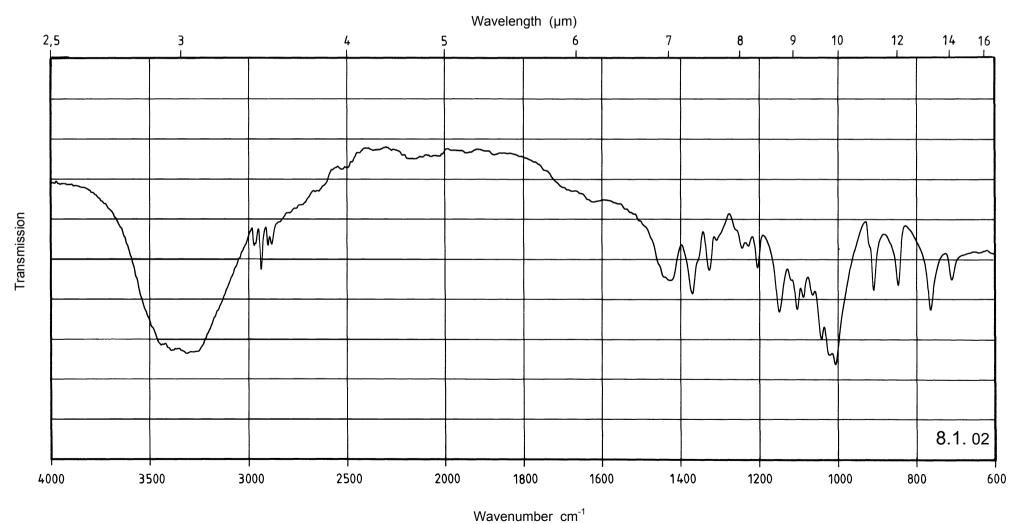


Name Fructose	Sample preparation Potassium bromide dispersion
	1.5 mg / 300 mg
	Reference
	Air
	Ordinate
	0 - 100 %

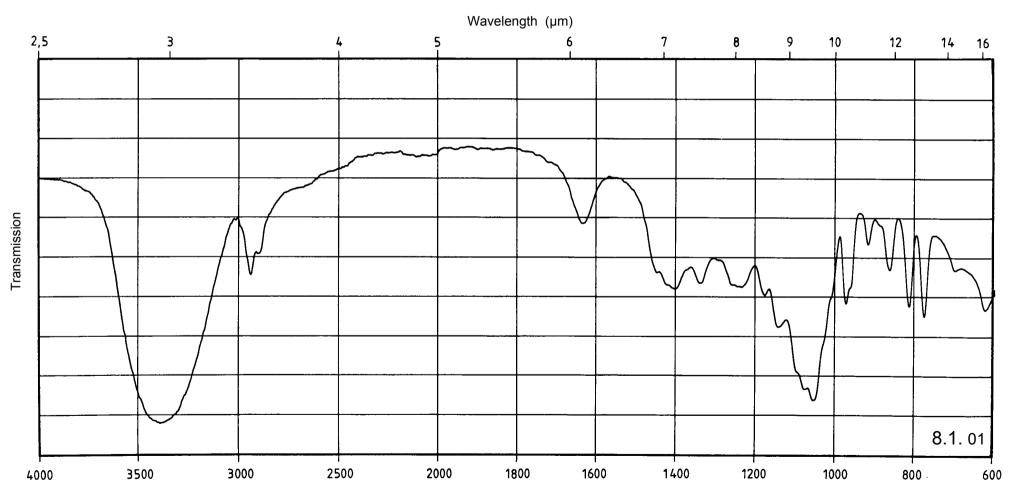


Name Lanolin, ethoxylated and acetylated Wool fat, ethoxylated and acetylated	Sample preparation Capillary film, sodium chloride cell
Solulan [®] 98	Reference Air Ordinate 0 - 100 %

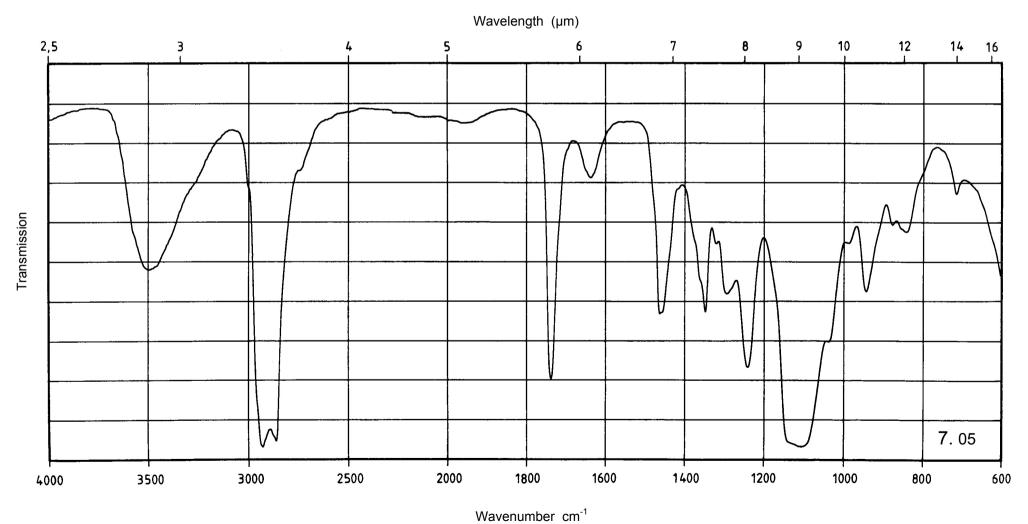


Name	Wool fat, Lanolin	Sample preparation Film, sodium chloride cell
		Reference
		Air
		Ordinate
		0 - 100 %

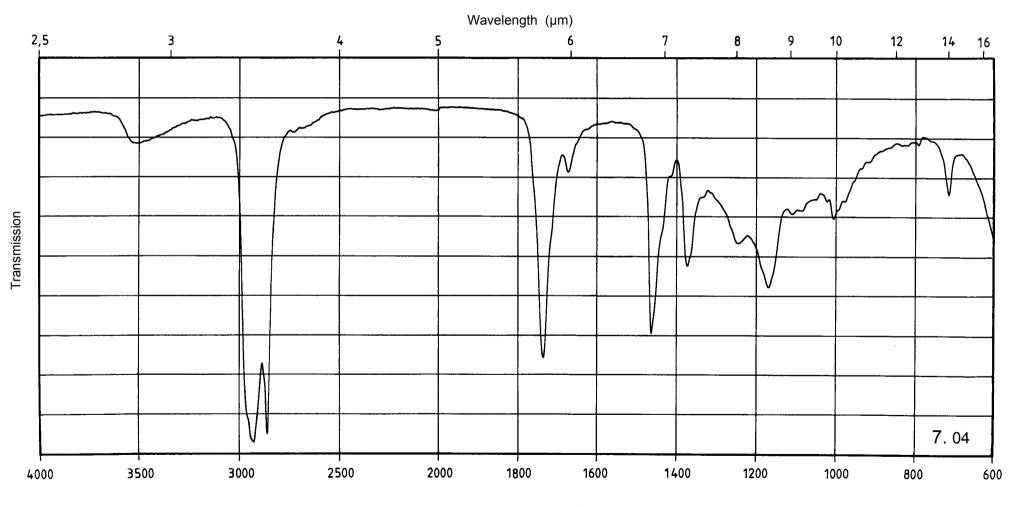


Name	Montan wax, ethanediol ester	Sample preparation
Ivanic		Potassium bromide dispersion
		3 mg / 300 mg
	Hoechst - Wachs E	Reference
	HOECHSL - Wachs E	Air
		Ordinate
		0 - 100 %

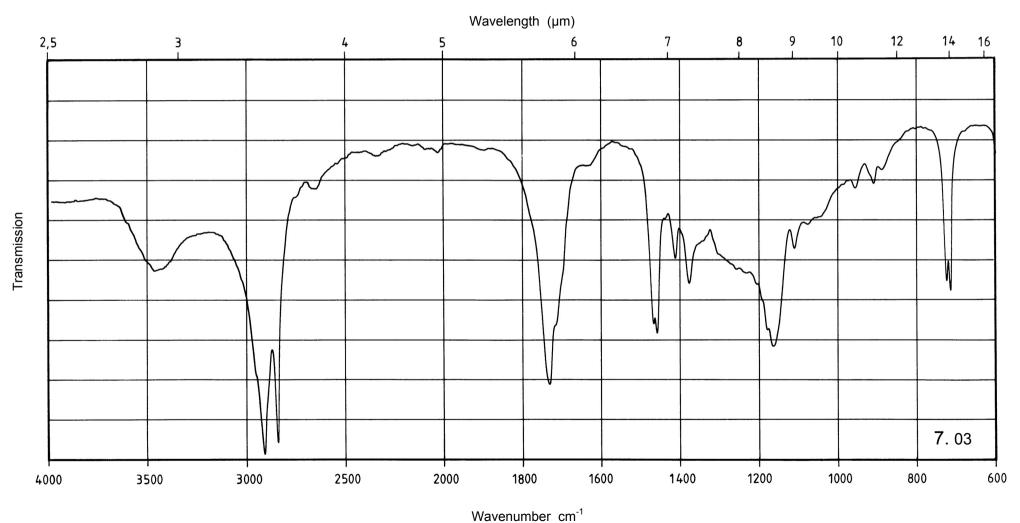


Name	Carnauba wax	Sample preparation Potassium bromide dispersion
	Brazil wax	2 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

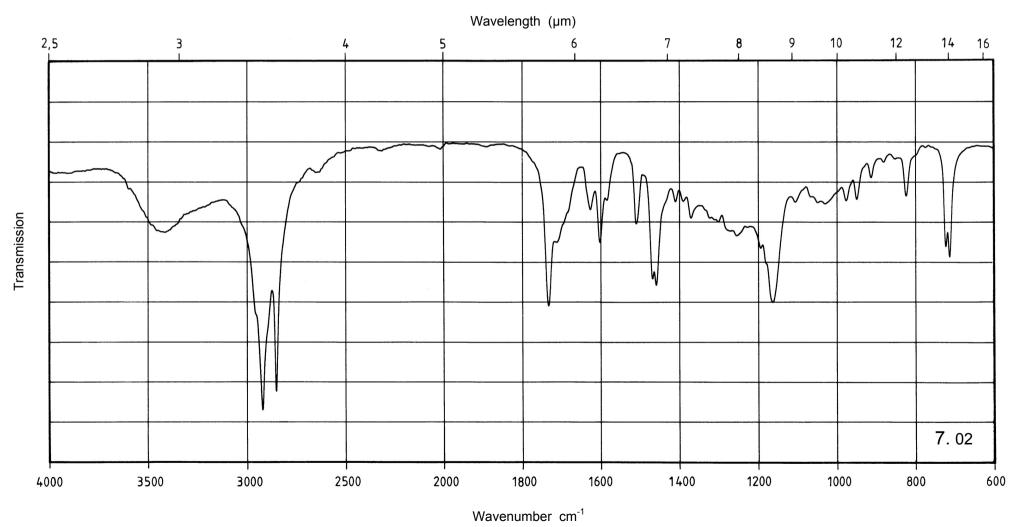


Name	Beeswax	Sample preparation Film, sodium chloride cell
		Reference
		Air
		Ordinate
		0 - 100 %

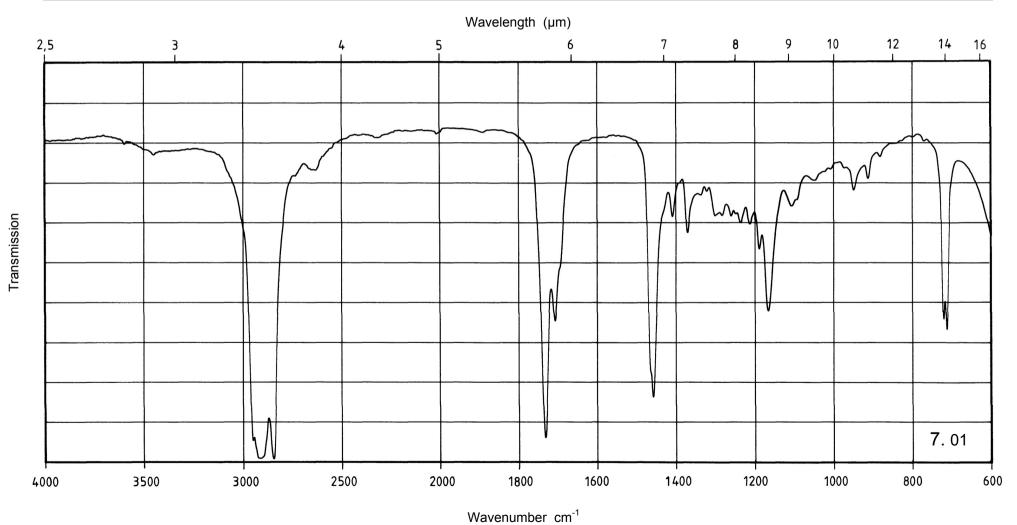


Name	Petrolatum	Sample preparation Film, sodium chloride cell
		Reference
		Air
		Ordinate
		0 - 100 %

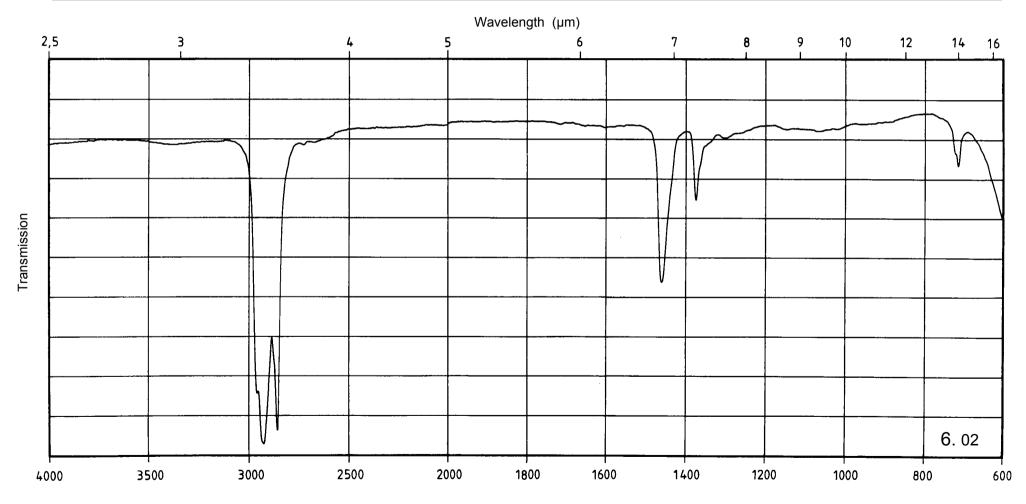


Name	Mineral Oil	Sample preparation Capillary film, sodium chloride cell
		Reference Air
		Ordinate 0 - 100 %

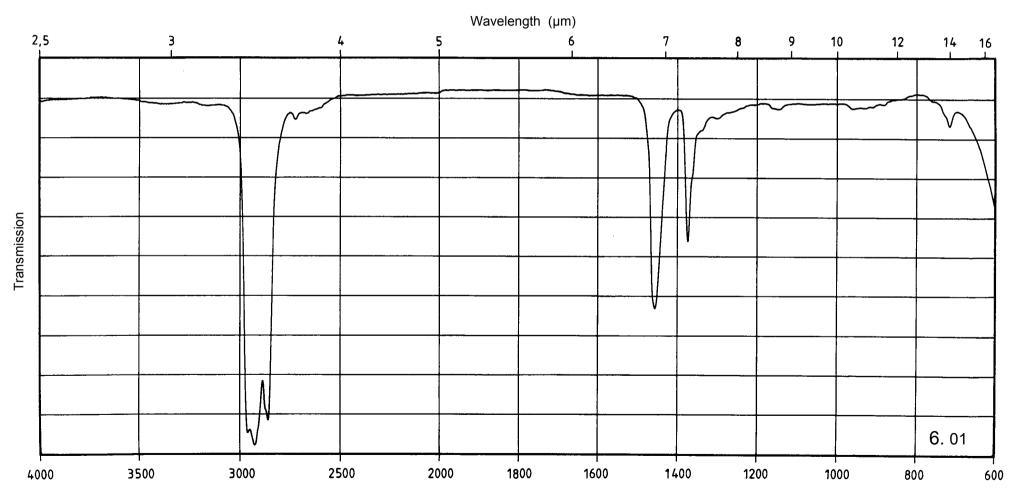


Name	Polyoxyethylene soy sterol	Sample preparation Capillary film, sodium chloride cell
	Generol® 122E10	Reference Air
		Ordinate 0 - 100 %

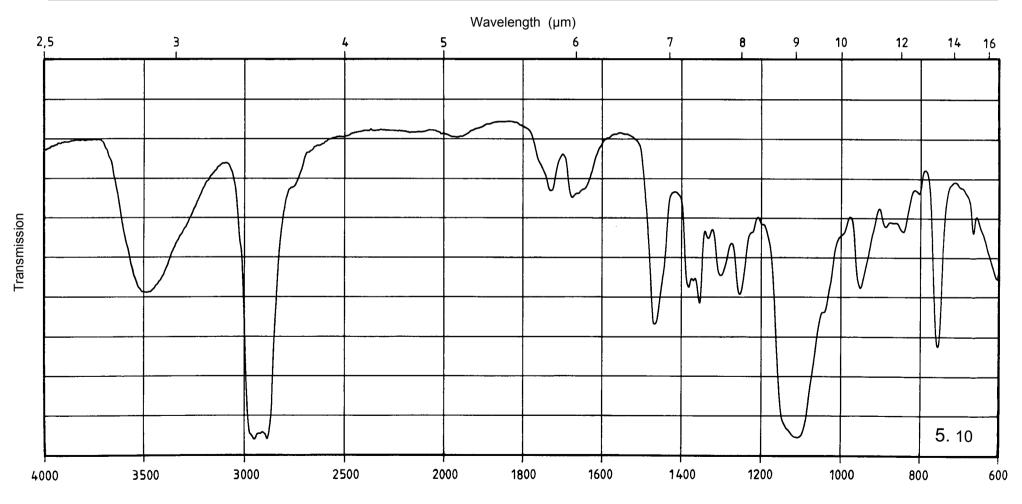


Name	Lanolin alcohol acetate	Sample preparation Capillary film, sodium chloride cell
	Acetulan ®	Reference Air
		Ordinate 0 - 100 %

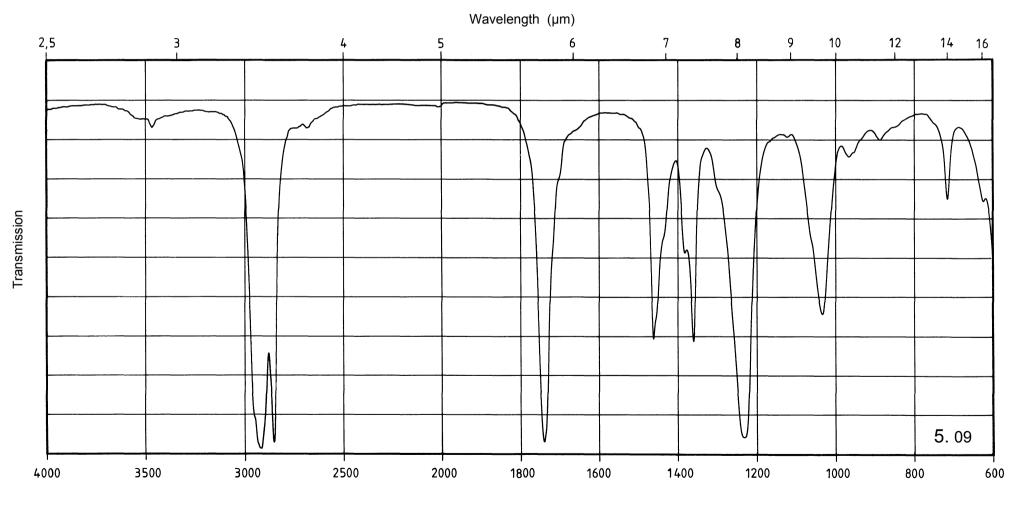


Name	Stearyl alcohol	Sample preparation
		Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

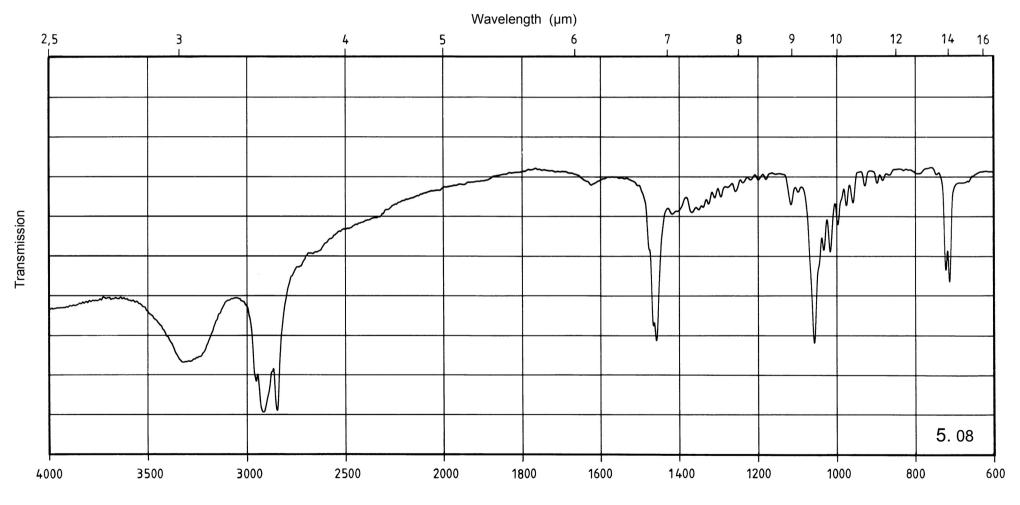


Name	Polyoxyethylene lauryl ether	Sample preparation Film, potassium bromide disk
	Brij [®] 35	Reference Air
		Ordinate 0 - 100 %

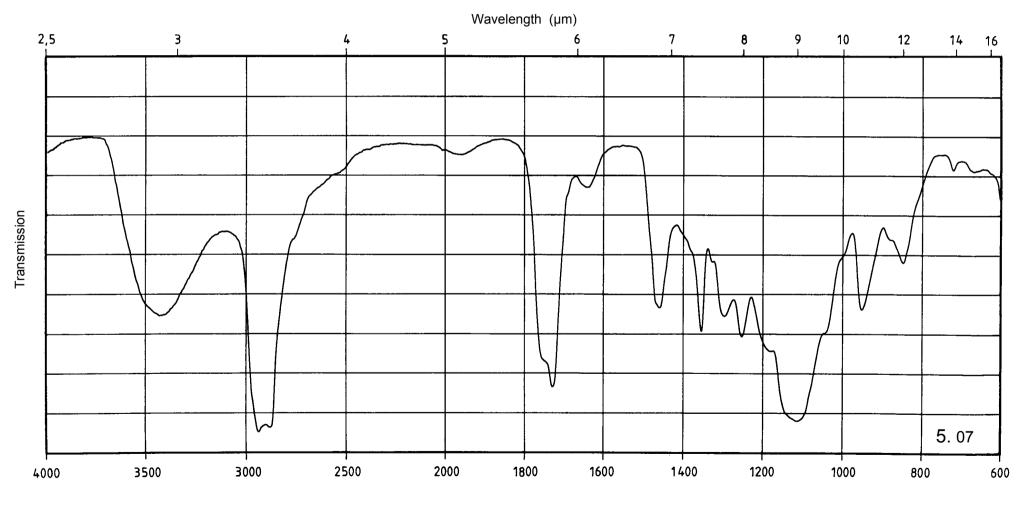


Name	2-Octyldodecanol	Sample preparation Capillary film, sodium chloride cell
	Eutanol® G	Reference Air
		Ordinate 0 - 100 %

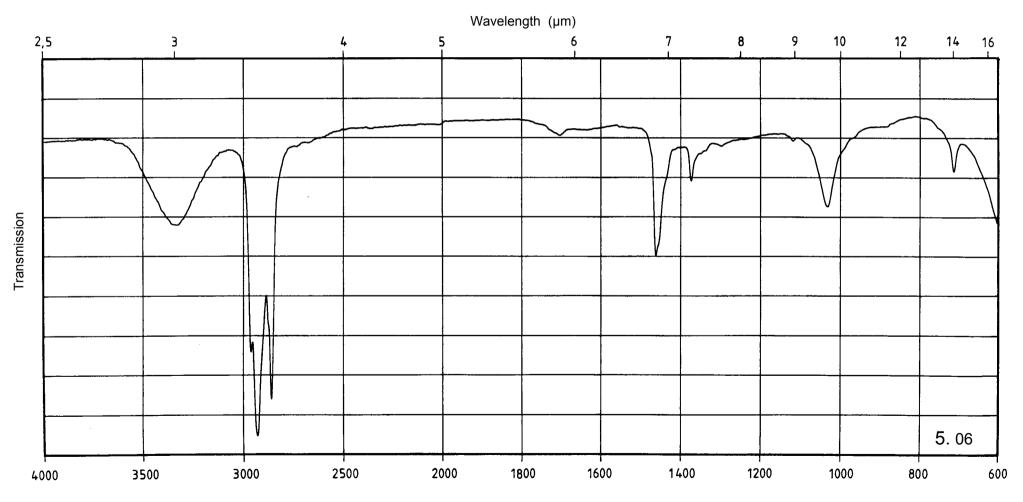


Name	Myristyl alcohol	Sample preparation Potassium bromide dispersion
		0.7 mg / 300 mg
	Lorol [®] C14	Reference Air
		Ordinate 0 - 100 %

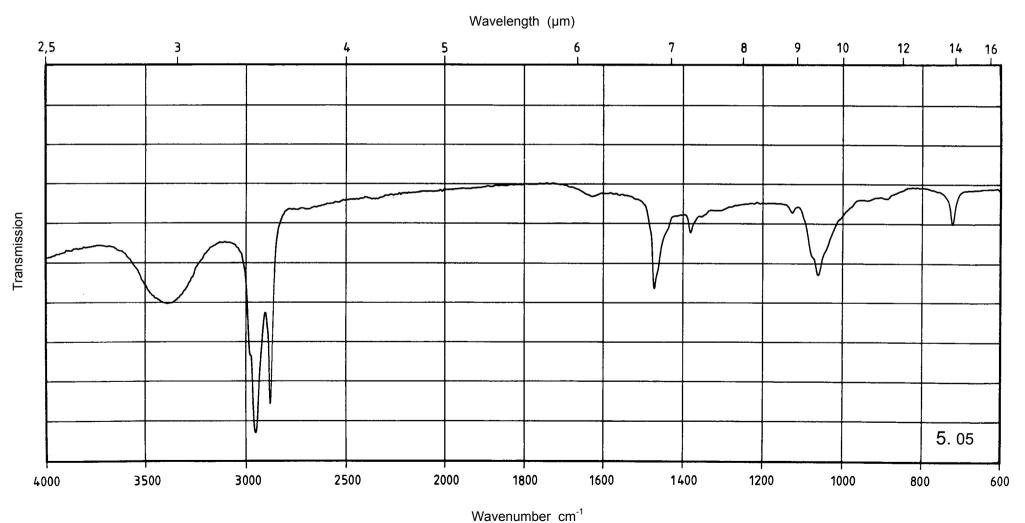


Name	Coconut, fatty alcohols	Sample preparation Film, sodium chloride cell
		Reference Air
		Ordinate 0 - 100 %

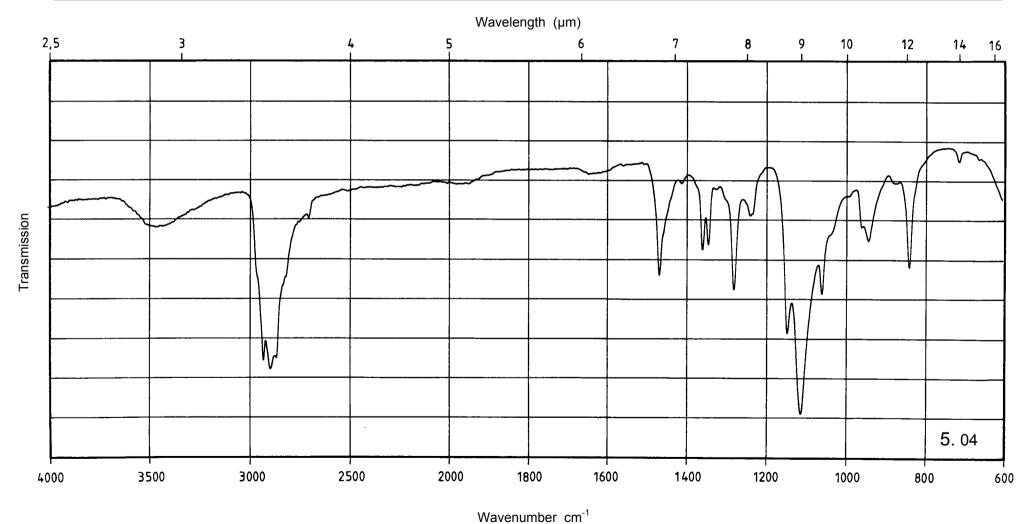


Name	Alcohols, fatty, ethoxylated	Sample preparation
		Potassium bromide dispersion
		1.5 mg / 300 mg
	Cremophor [®] O	Reference
	Cremophor O	Air
		Ordinate
		0 - 100 %

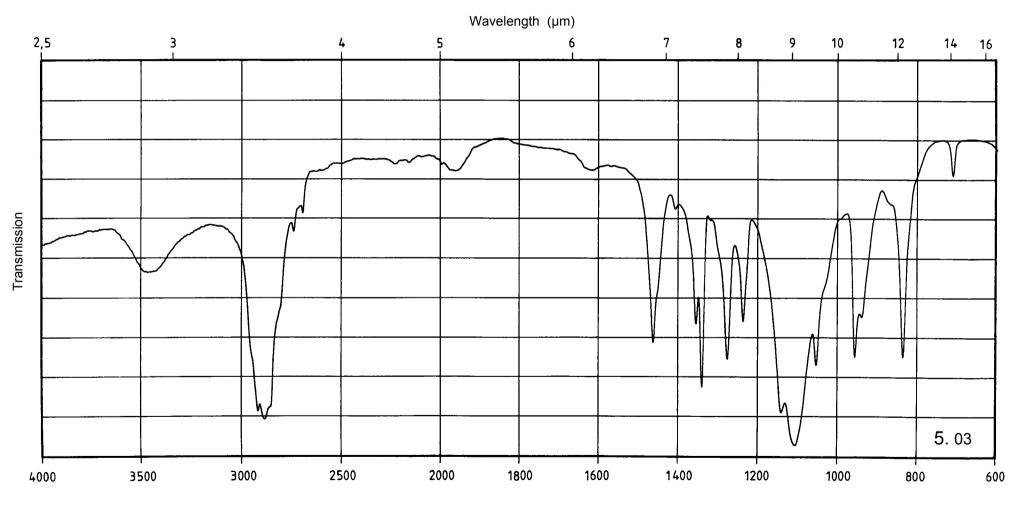


Name Cetyl stearyl alco	Cetyl stearyl alcohol, sodium cetyl stearyl sulfate, 9 + 1	Sample preparation
Ivanic	deligible octyl stearyr alcorror, socialli cetyl stearyr sallate, 5 i i	Potassium bromide dispersion
		1 mg / 300 mg
	Lanette ® N	Reference
	Lanette in	Air
		Ordinate
		0 - 100 %

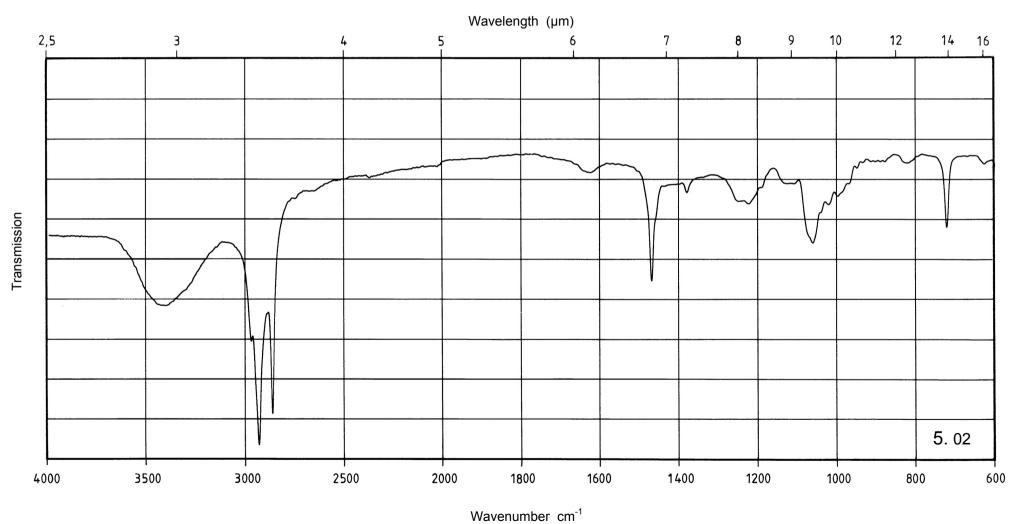


Name	Cetyl alcohol	Sample preparation
		Potassium bromide dispersion
		1 mg / 300 mg
	Lorol® C16	Reference
		Air
		Ordinate
		0 - 100 %

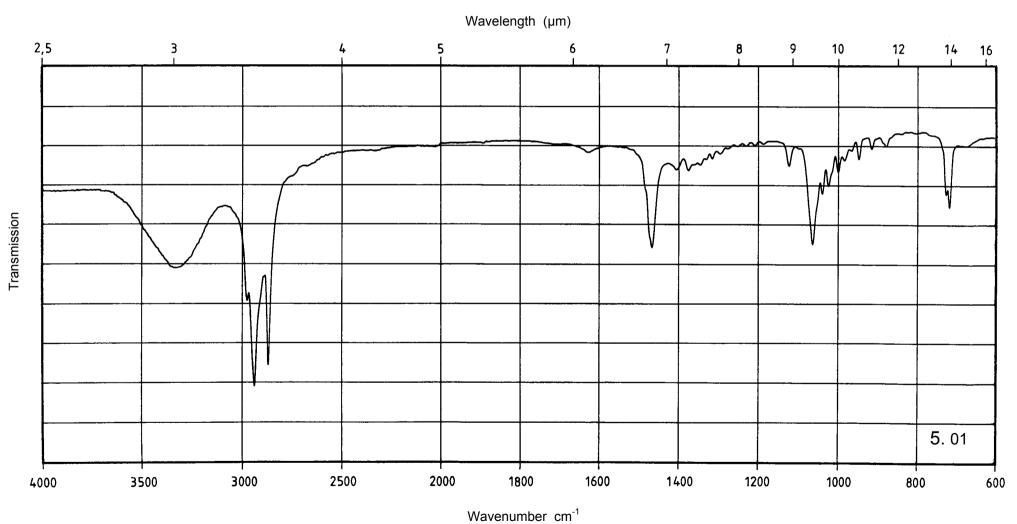


Name	12-Hydroxystearic acid polyglycol ester	Sample preparation Capillary film, sodium chloride cell
	Solutol [®] HS 15	Reference Air
		Ordinate 0 - 100 %

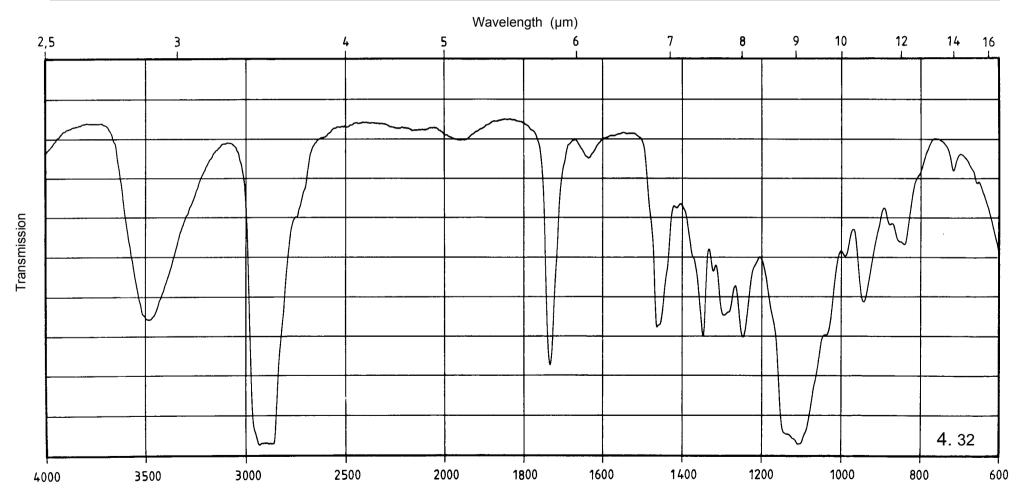


Name	Triglycerides (C ₈ - C ₁₂)	Sample preparation Capillary film, sodium chloride cell
	Miglyol® 812	Reference Air
		Ordinate 0 - 100 %

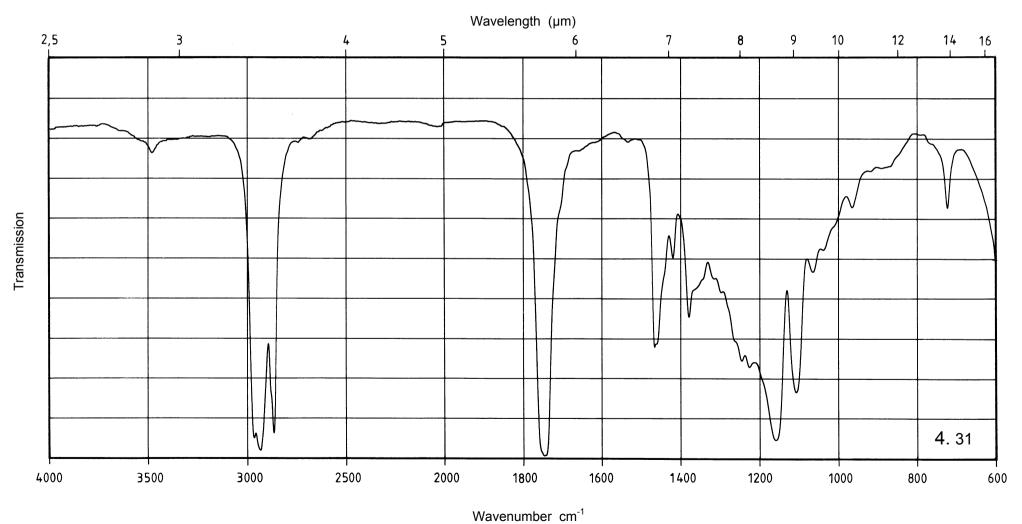


Name	Tegacid [®] H, Tegenacid [®] H	Sample preparation Potassium bromide dispersion
		0.8 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

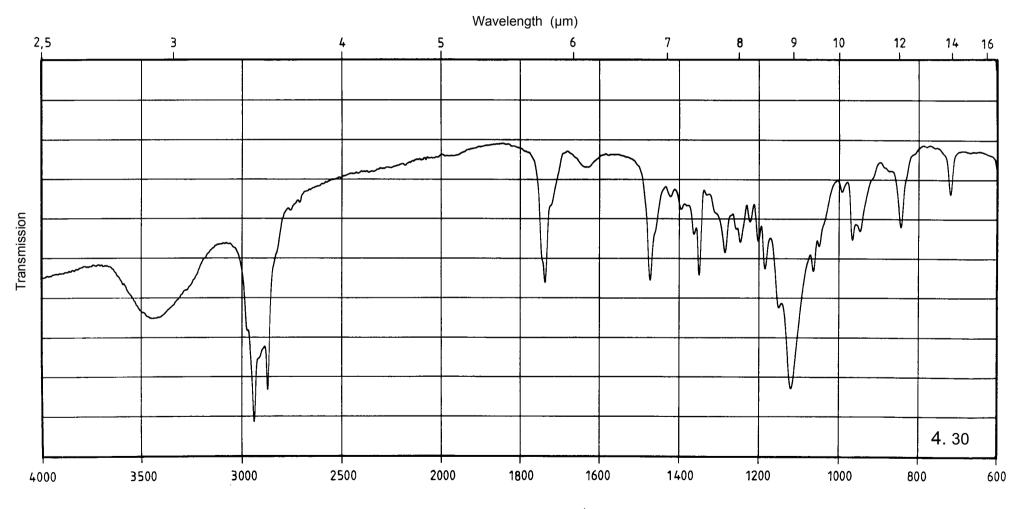


Name	Caprylic/capric/stearic triglyceride	Sample preparation Capillary film, sodium chloride cell
	Softisan® 378	Reference Air
		Ordinate 0 - 100 %

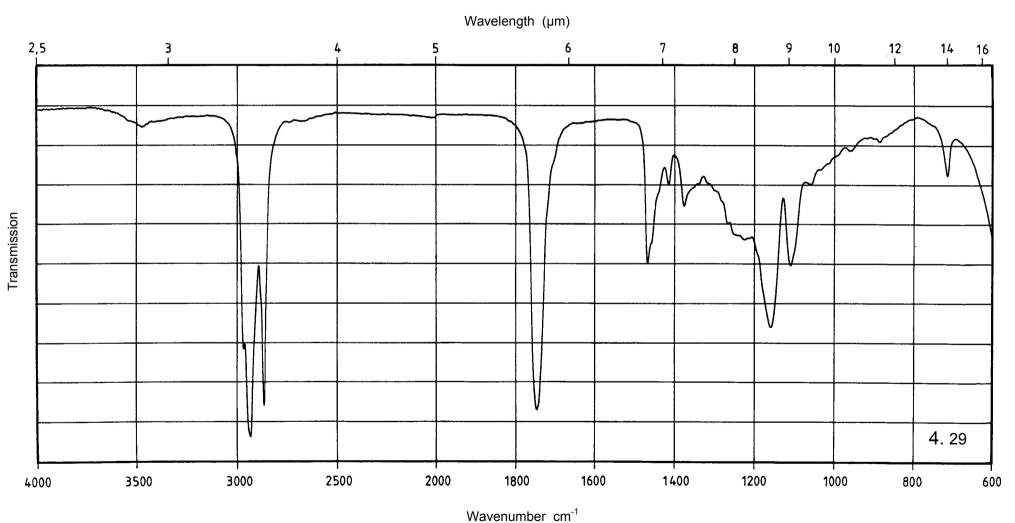



Name	Stearic acid	Sample preparation
		Potassium bromide dispersion
		1 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

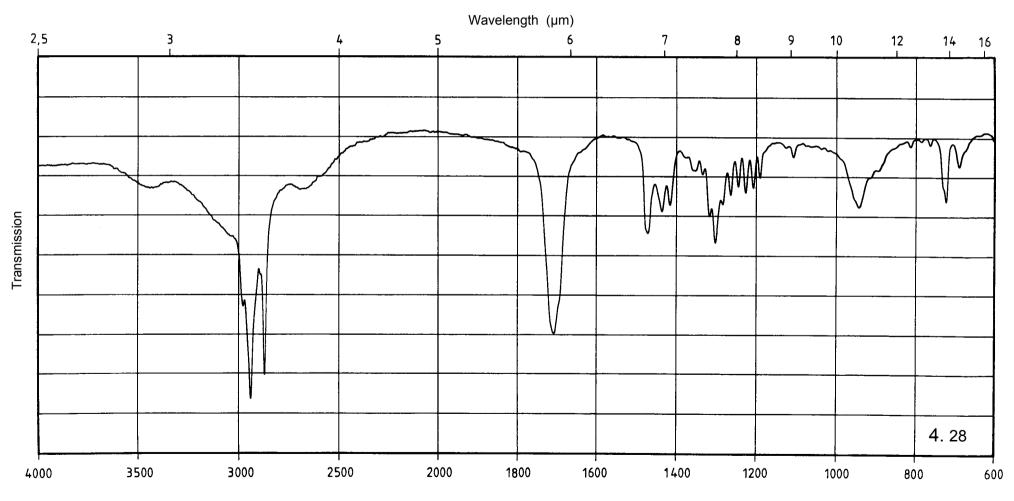


Name	Stearin	Sample preparation
Ivanic		Potassium bromide dispersion
		1 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

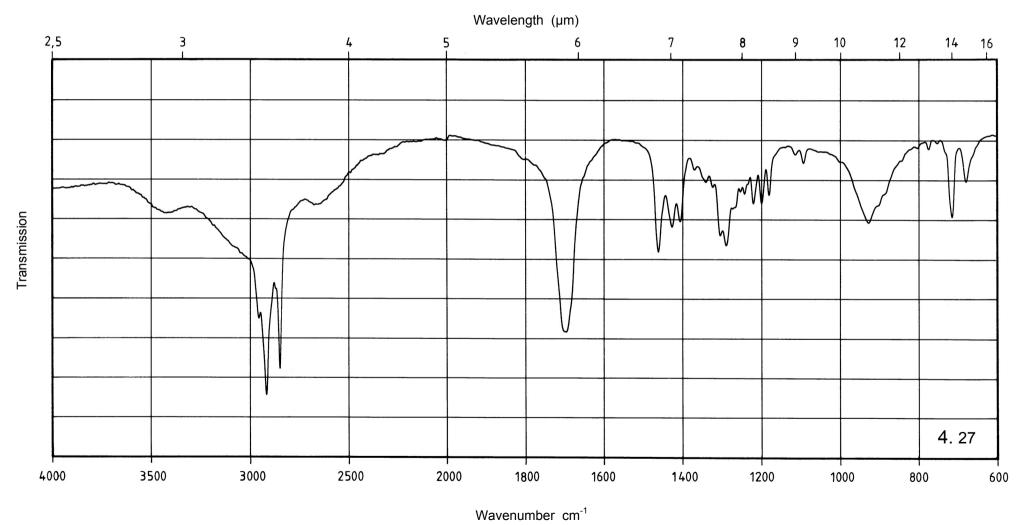

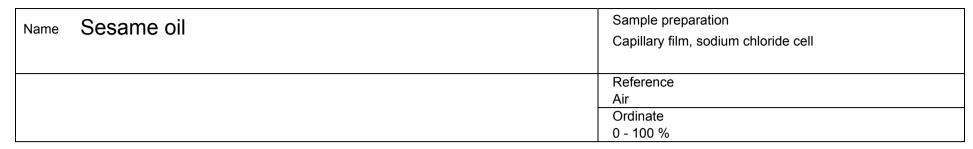


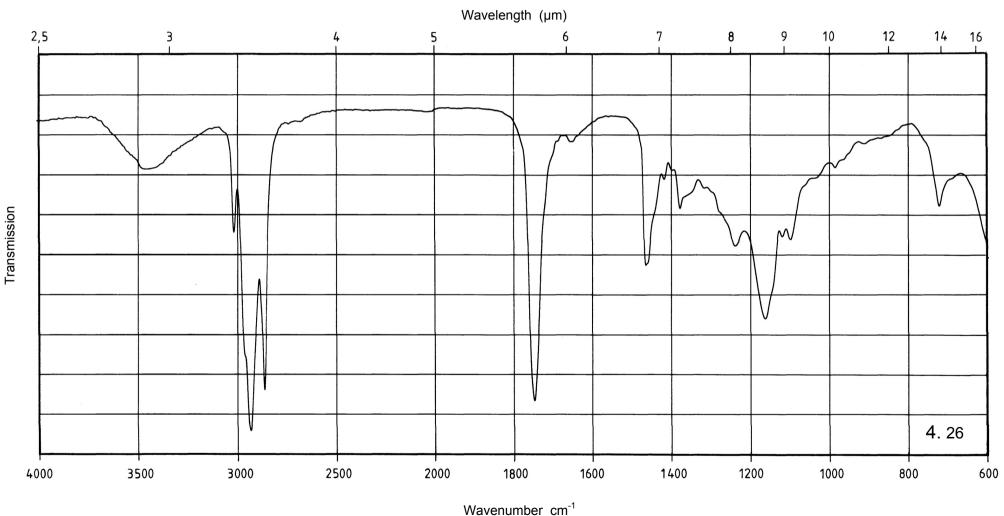


Name	Castor oil, hydrogenated, ethoxylated	Sample preparation Capillary film, sodium chloride cell
	Cremophor® RH 60	Reference Air
		Ordinate 0 - 100 %

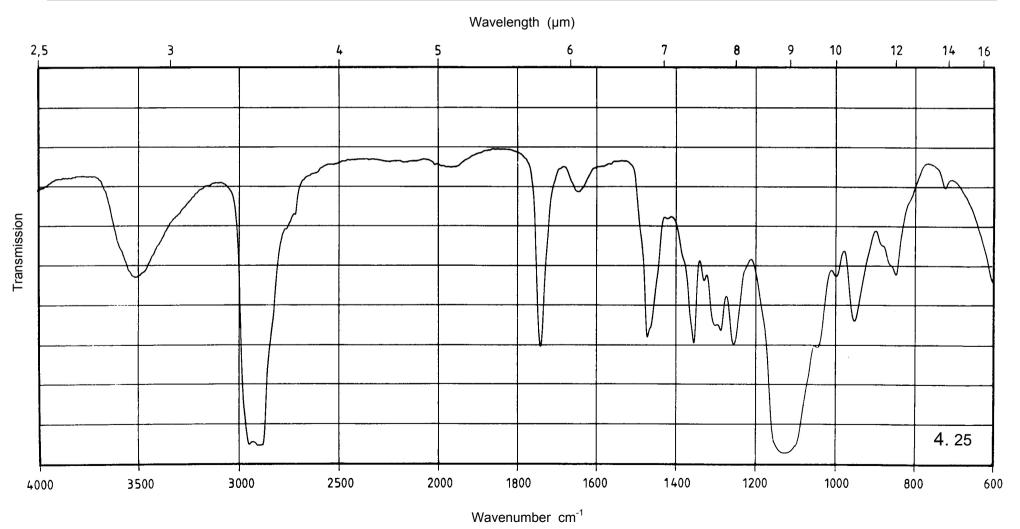


Name	Castor oil, ethoxylated	Sample preparation Capillary film, sodium chloride cell
	Cremophor [®] EL	Reference Air
		Ordinate 0 - 100 %

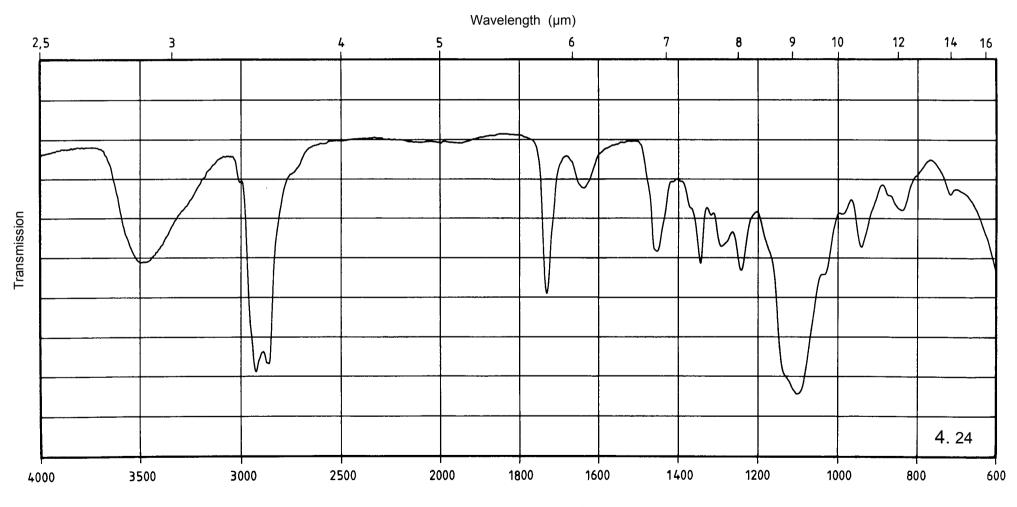


Sample preparation Capillary film, sodium chloride cell
Reference Air
Ordinate 0 - 100 %

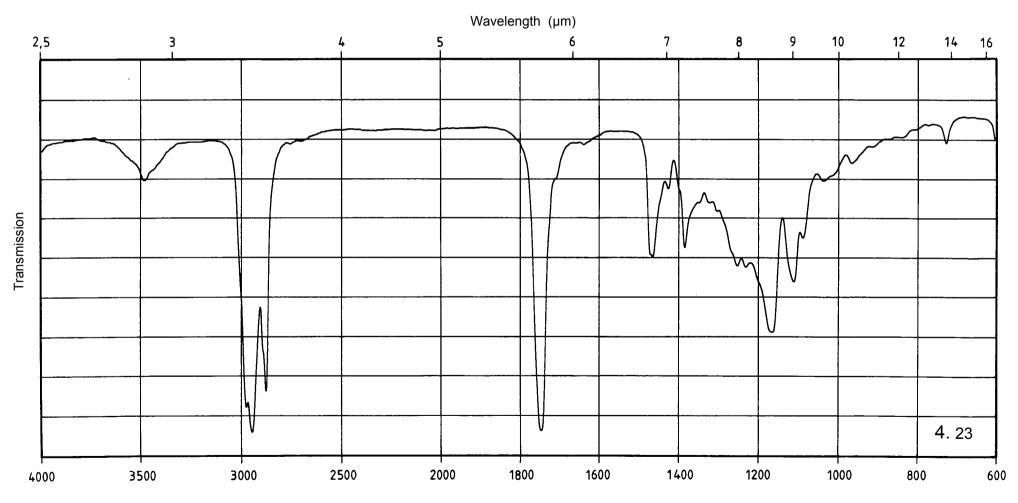


Name	Polyoxyethylene 50 stearate	Sample preparation
Ivanic		Potassium bromide dispersion
		1 mg / 300 mg
	Myrj [®] 53	Reference
	Myrj 55	Air
		Ordinate
		0 - 100 %

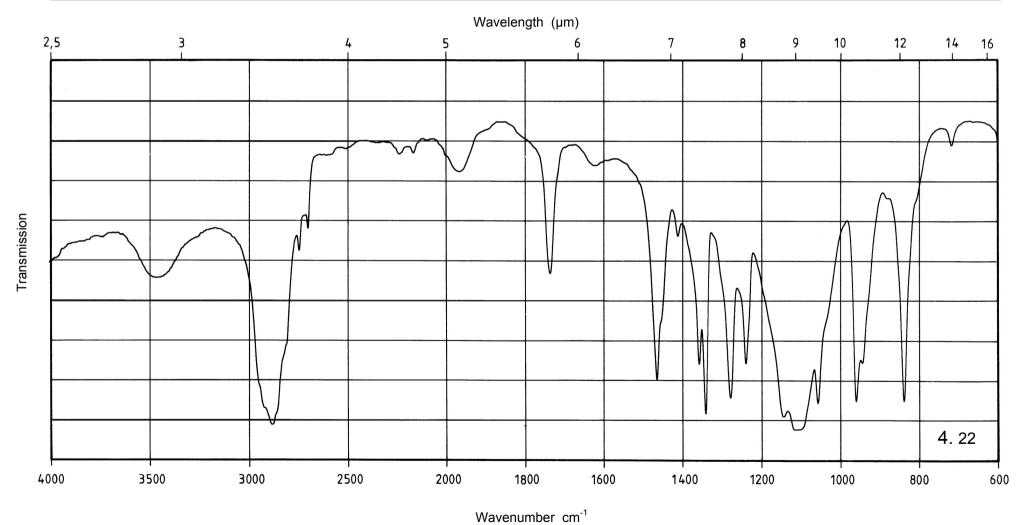


Name	Polyoxyethylene fatty acid ester	Sample preparation Capillary film, sodium chloride cell
	Arlacel® 989	Reference Air
		Ordinate 0 - 100 %

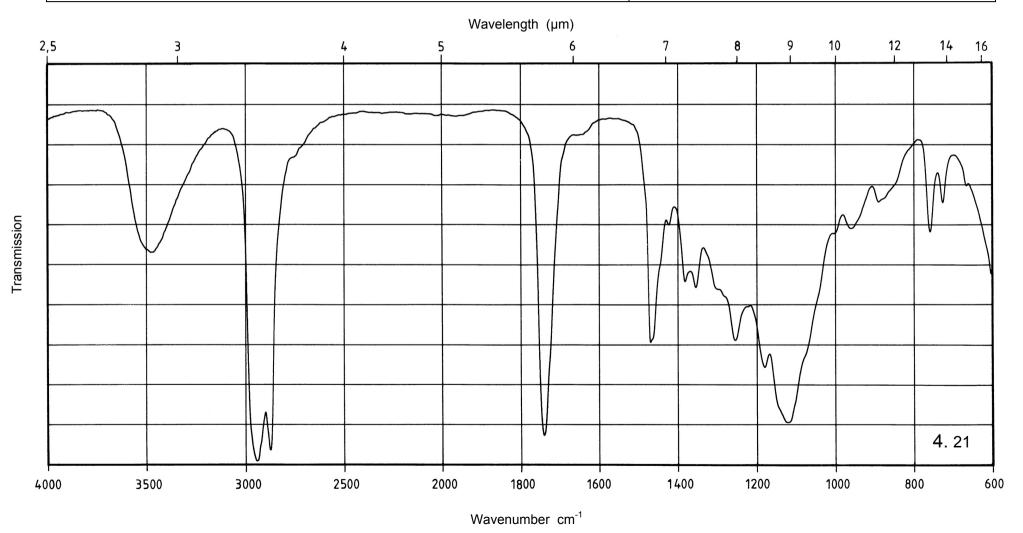


Name	Poly(oxyethylene) glycerol fatty acid ester	Sample preparation Capillary film, sodium chloride cell
	Cetiol® HE	Reference Air
		Ordinate 0 - 100 %

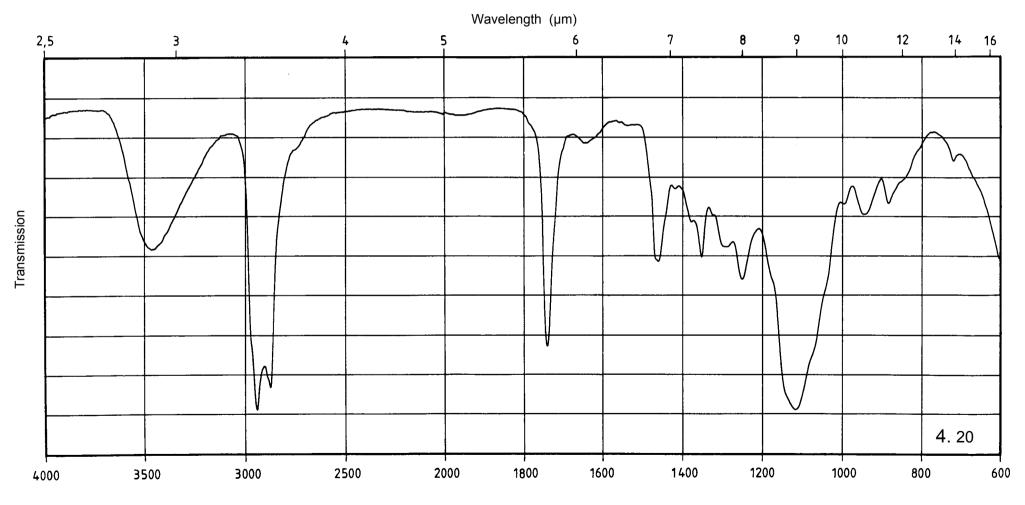


Name	Palmitic acid	Sample preparation
Name		Potassium bromide dispersion
	Hexadecanoic acid	1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

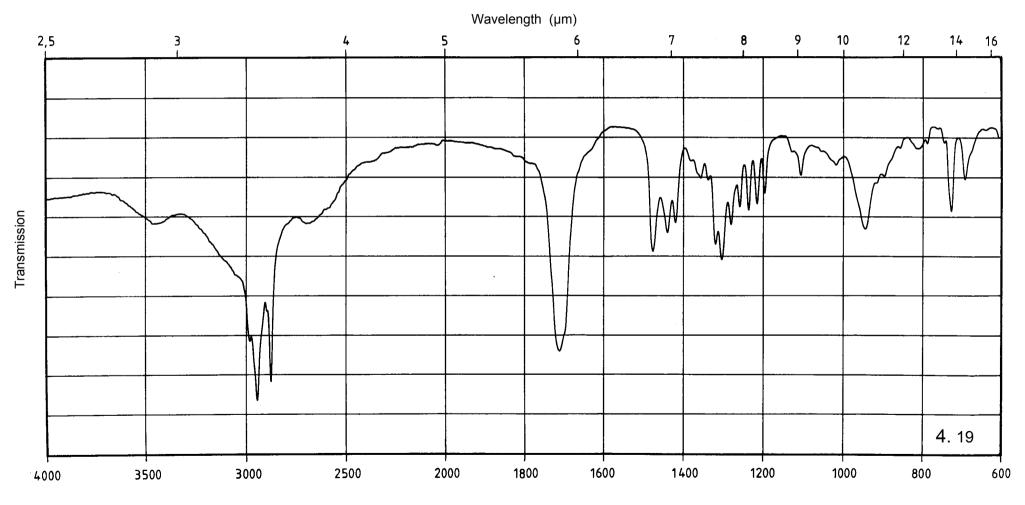


Name	Oleic acid, decyl ester	Sample preparation Capillary film, sodium chloride cell
	Cetiol® V	Reference Air
		Ordinate 0 - 100 %

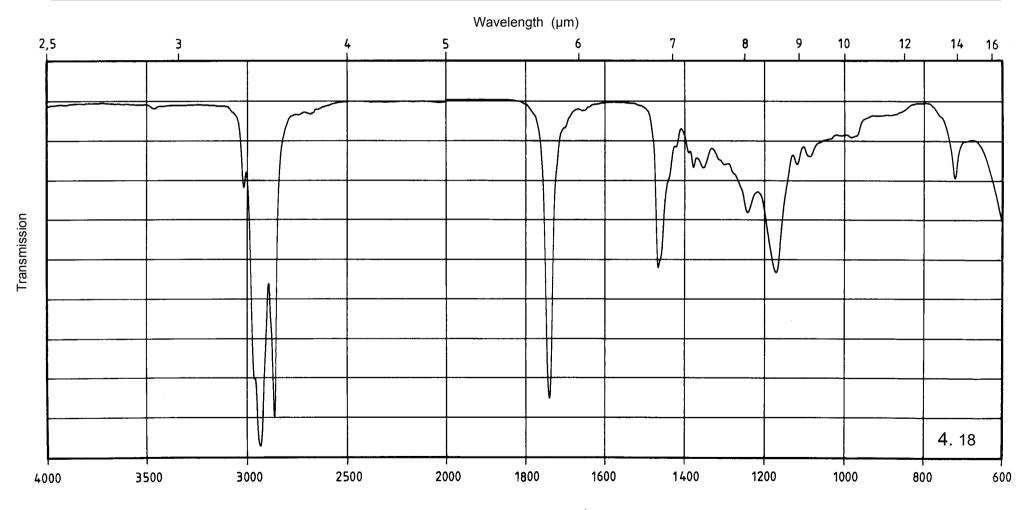



Name	Octyl palmitate	Sample preparation Capillary film, sodium chloride cell
	Ceraphyl® 368	Reference Air
		Ordinate
		0 - 100 %

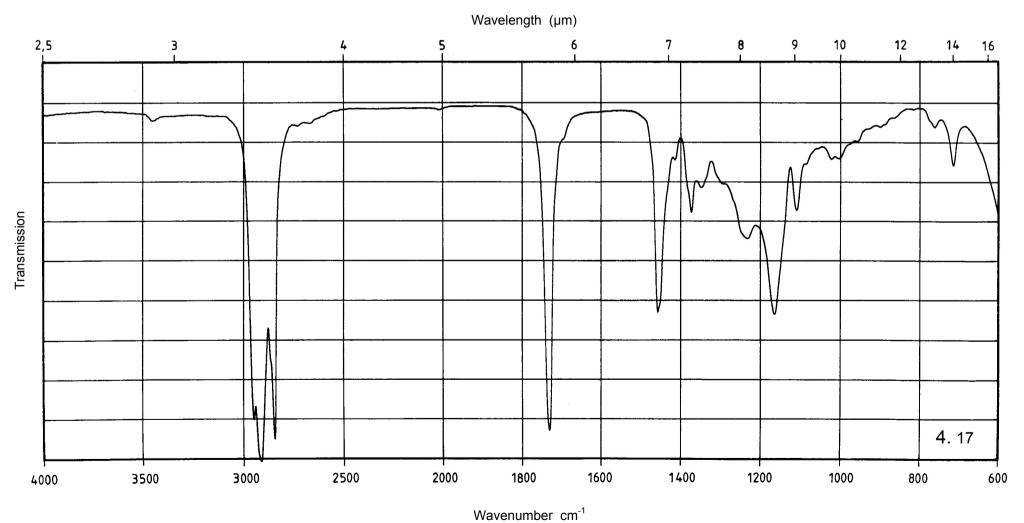

Name	Sodium stearyl fumarate	Sample preparation
		Potassium bromide dispersion
		0.9 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

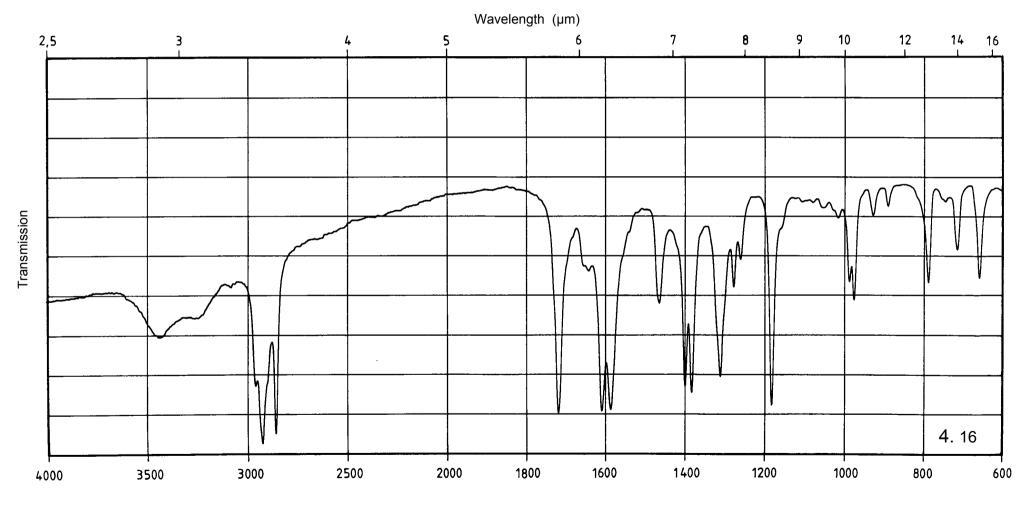


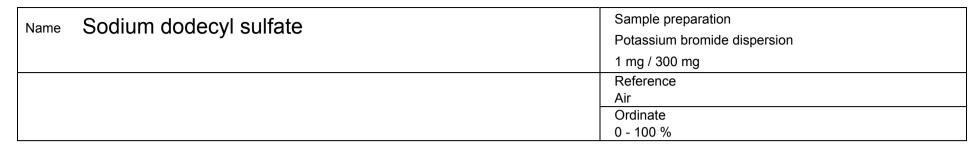

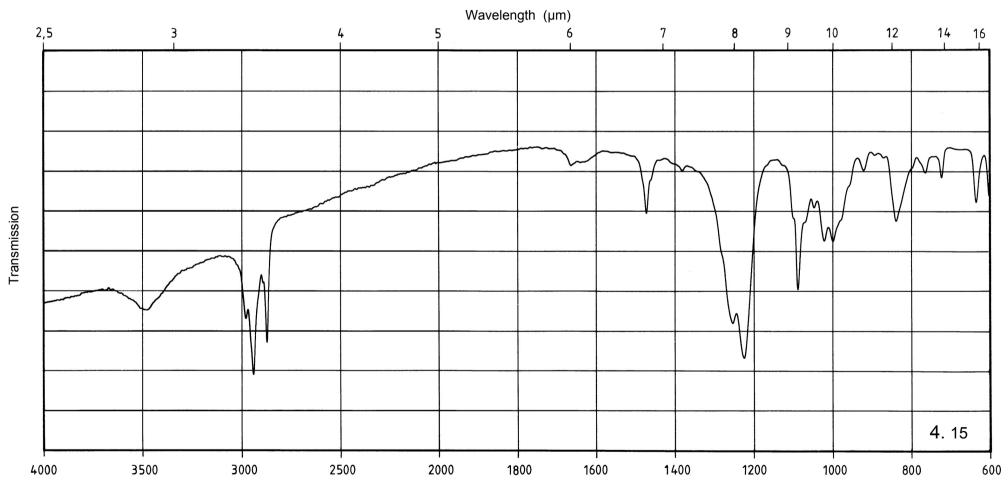


Wavenumber cm⁻¹


Name	Monoglycerides, acetylated	Sample preparation Capillary film, sodium chloride cell
	Cetodan® 90 - 40, Myvacet®	Reference Air
		Ordinate 0 - 100 %

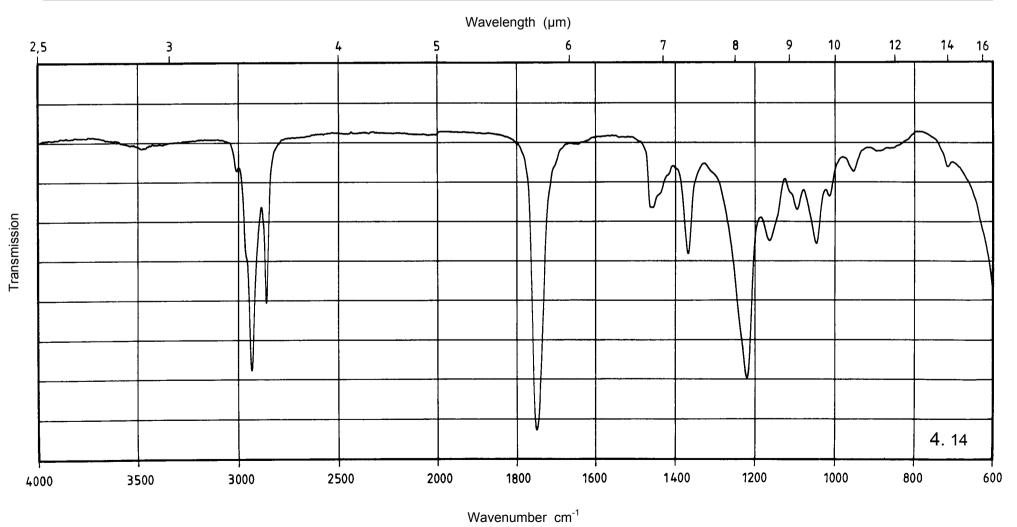


Name	Methyl glucose sesquistearate	Sample preparation
Name Wictify! gladose sesquistedrate	Wetry glacose sesquistearate	Potassium bromide dispersion
		1.5 mg / 300 mg
	Glucate ® SS	Reference
	Glucale 33	Air
		Ordinate
		10 - 100 %

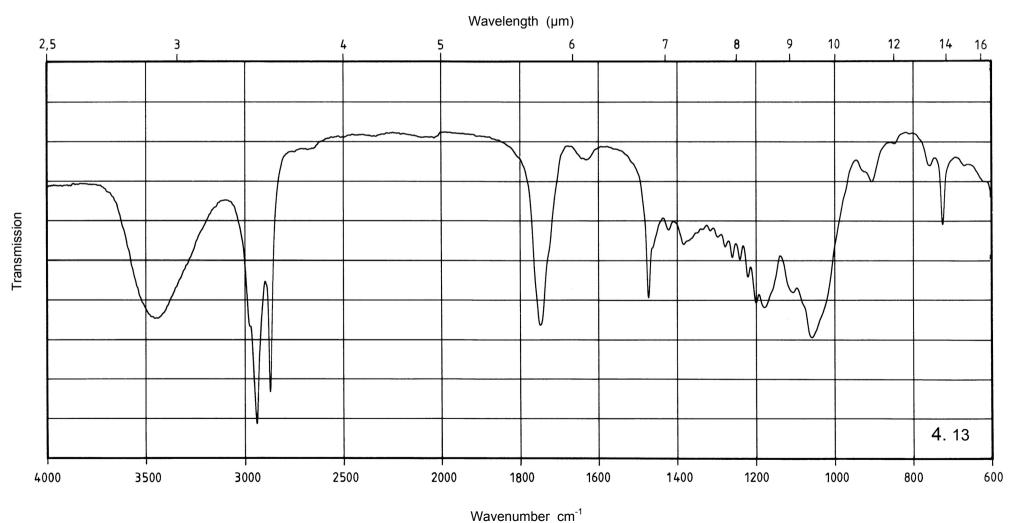


Name	Magnesium stearate	Sample preparation
		Potassium bromide dispersion
		0.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 80 %

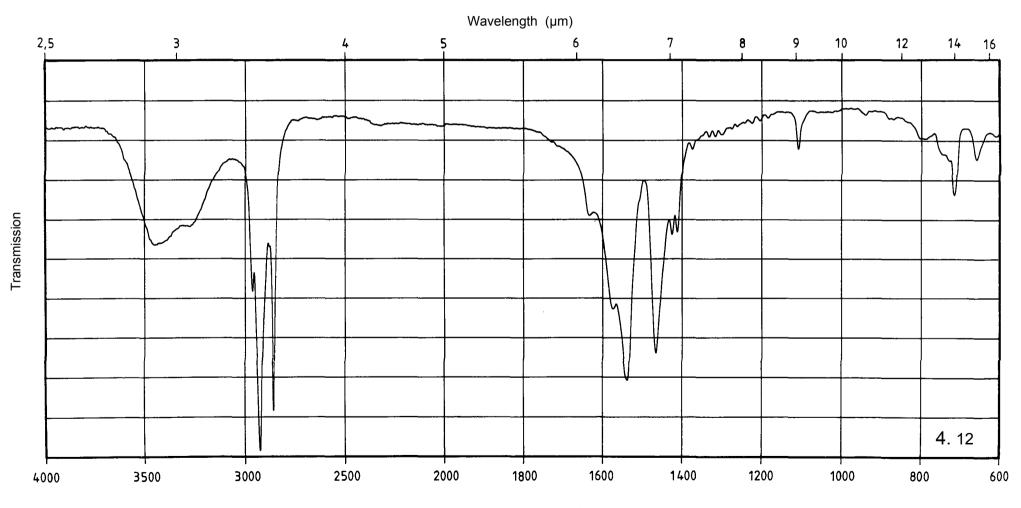


Name	Name Coconditiativ acid diethanolamide	Sample preparation
Ivanic		Potassium bromide dispersion
		1 mg / 300 mg
	Comperlan [®] KD	Reference
	Compenan KD	Air
		Ordinate
		10 - 90 %

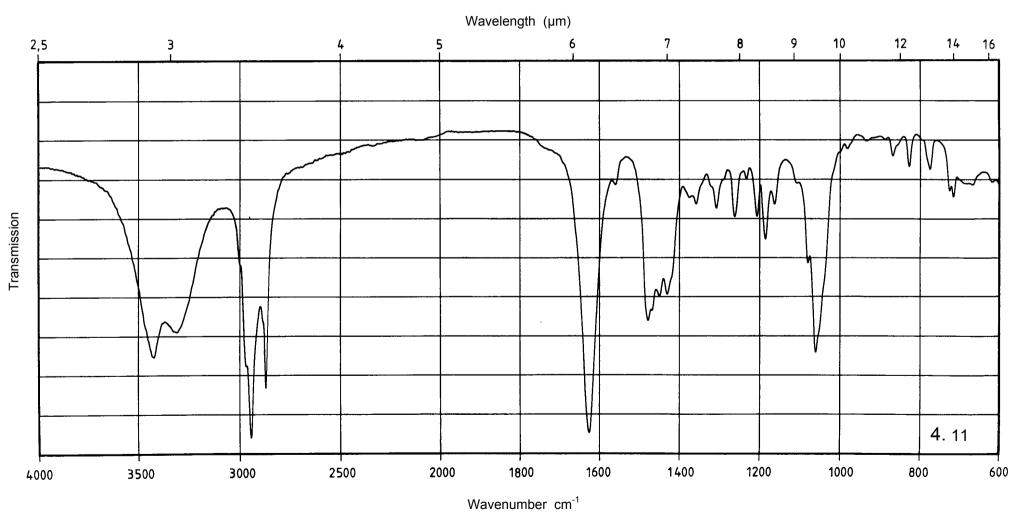


Name	Coconut oil	Sample preparation Capillary film, sodium chloride cell
		Reference Air
		Ordinate 0 - 100 %

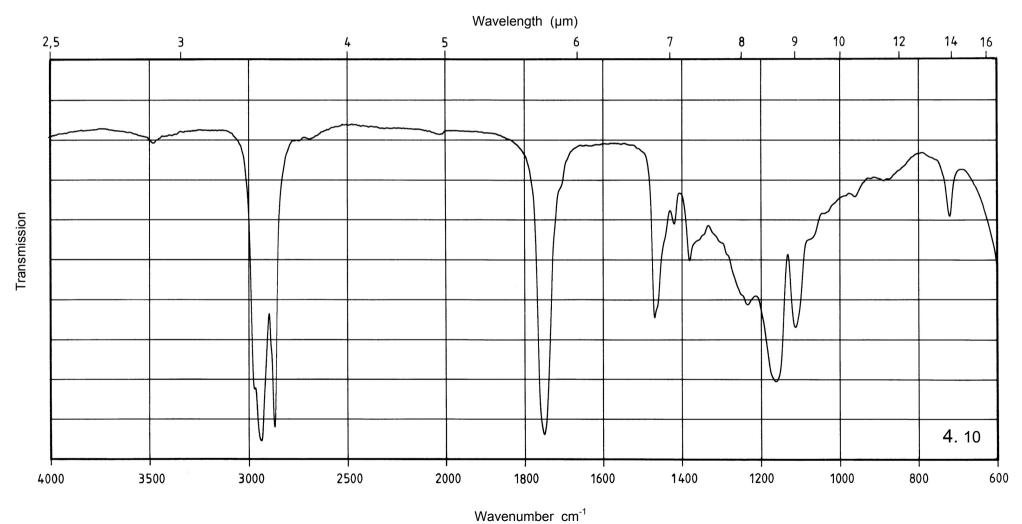


Name	11-Hydroxystearic acid	Sample preparation
Ivanic		Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

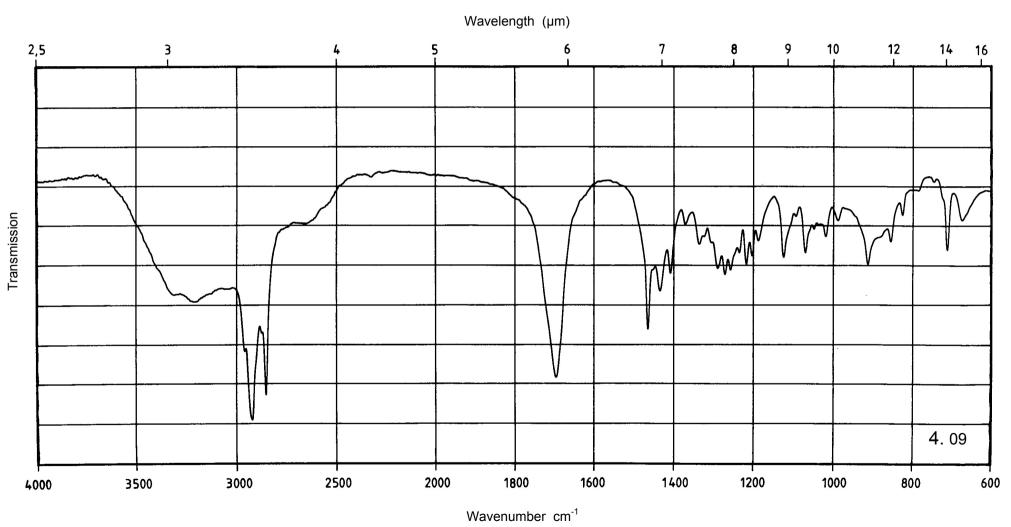


Name	Triglycerides (C ₁₂ - C ₁₈)	Sample preparation Film, potassium bromide disk
		Reference Air
		Ordinate 0 - 100 %

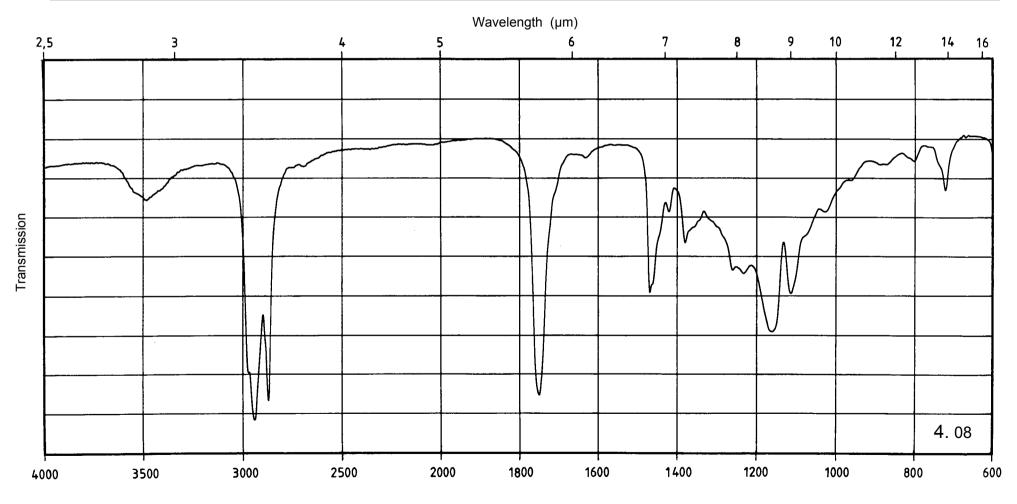


Name	Glycerol monostearate	Sample preparation
Name		Potassium bromide dispersion
		1 mg / 300 mg
	Cutina® GMS, Tegin® M	Reference
	Cullia Givio, regin ivi	Air
		Ordinate
		0 - 100 %

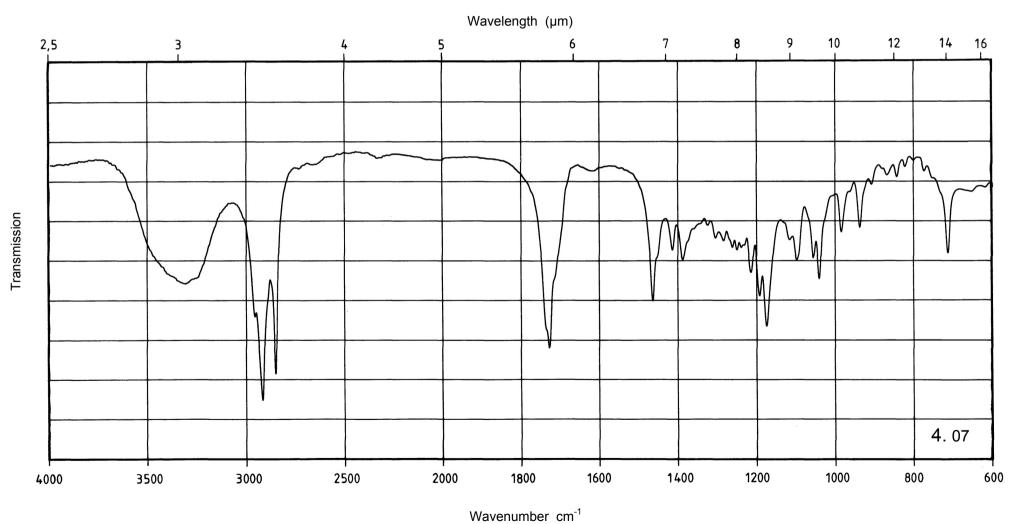


Name		Sample preparation
Ivanic		Potassium bromide dispersion
		1 mg / 300 mg
	Arlacel [®] 165	Reference
	Allacei 100	Air
		Ordinate
		0 - 100 %

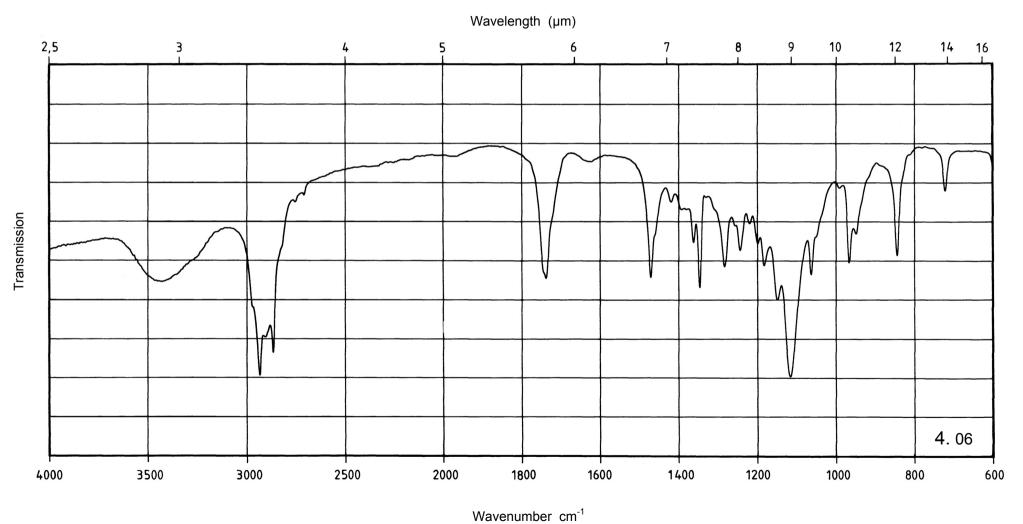



Name	Fatty acid sarcosides	Sample preparation Capillary film, sodium chloride cell
	Medialan [®] KF	Reference Air
		Ordinate 0 - 100 %

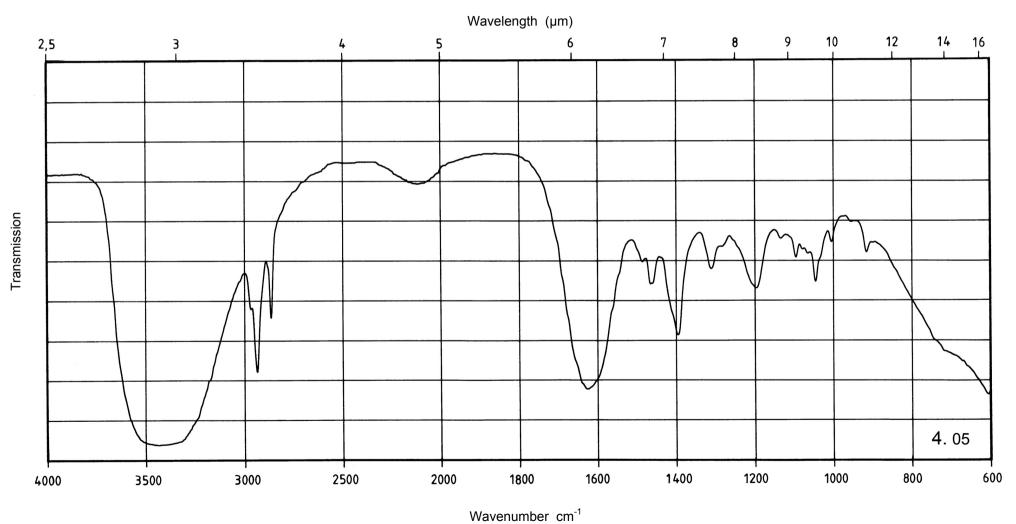


Name	Triglycerides, Mixture of	Sample preparation
		Potassium bromide dispersion
		1 mg / 300 mg
	Emulsan [®] MD	Reference
	ETTUSATI IVID	Air
		Ordinate
		0 - 100 %

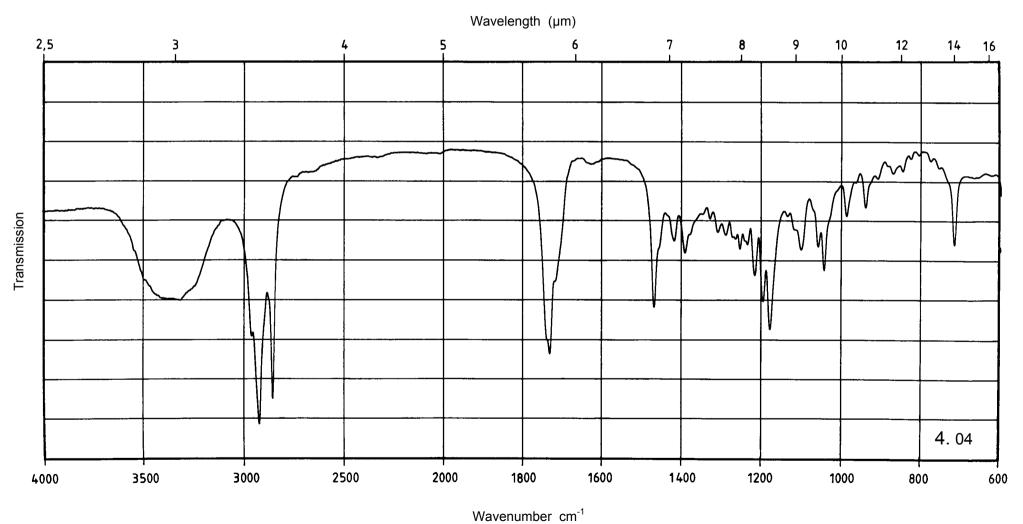

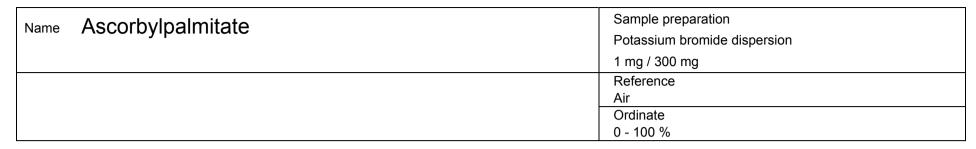


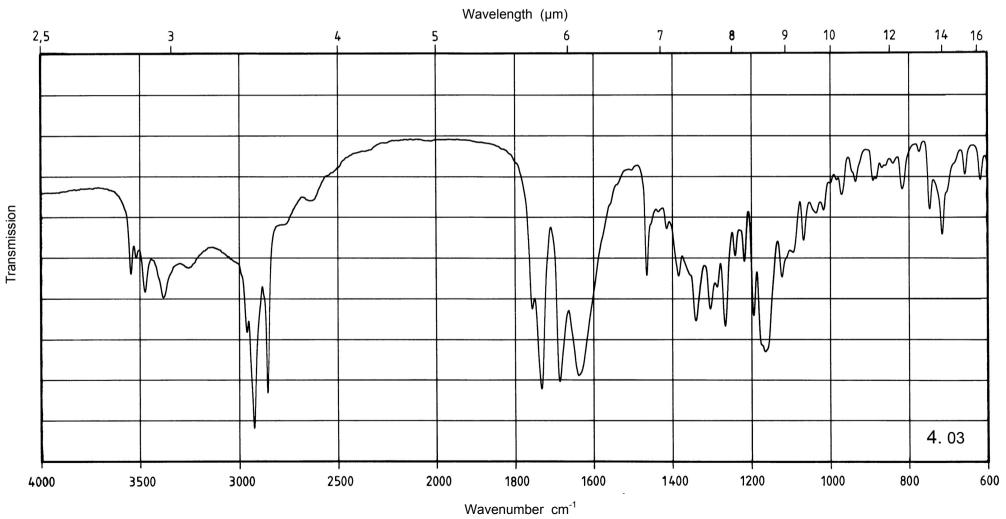


Name	Aluminium stearate	Sample preparation
		Potassium bromide dispersion
		1 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

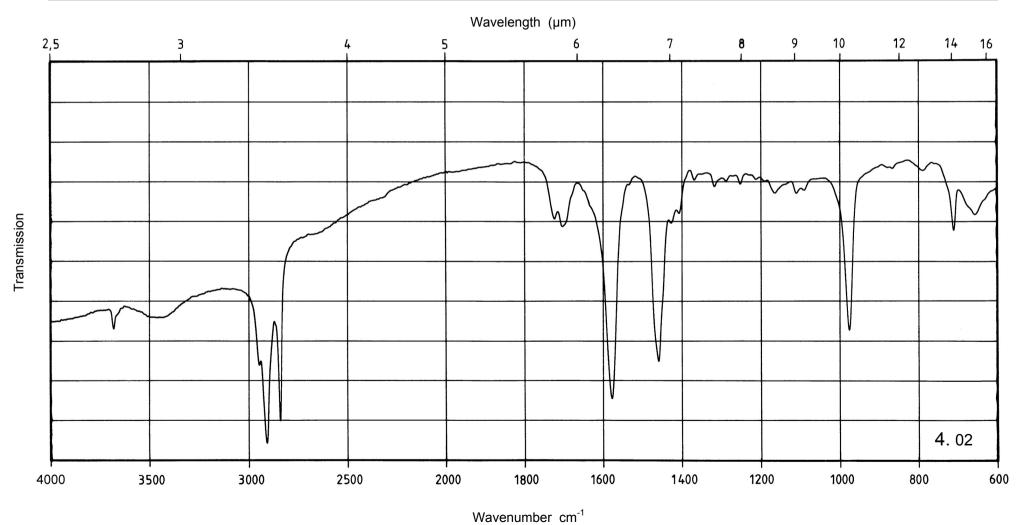


Name	Aluminium monostearate	Sample preparation
Ivanic		Potassium bromide dispersion
		1 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

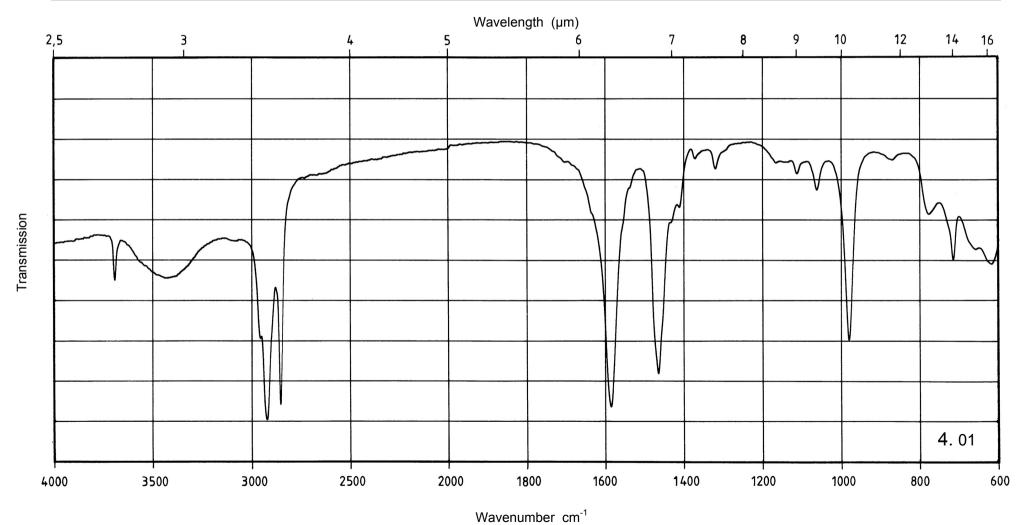


Name	Triethylcitrate	Sample preparation Capillary film, sodium chloride cell
	Citroflex [®]	Reference Air
		Ordinate 0 - 100 %

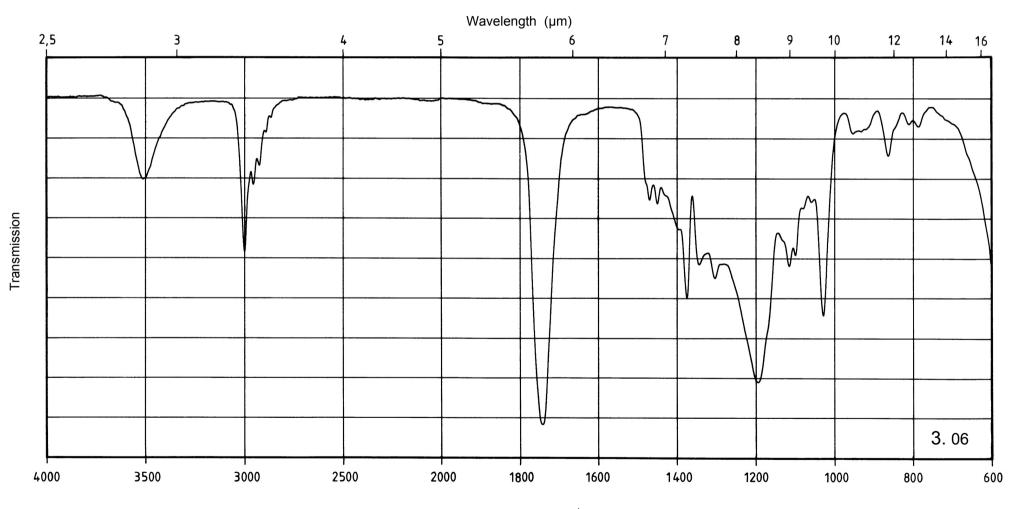


Name Triacetine	Sample preparation Capillary film, sodium chloride cell
	Reference Air
	Ordinate 0 - 100 %

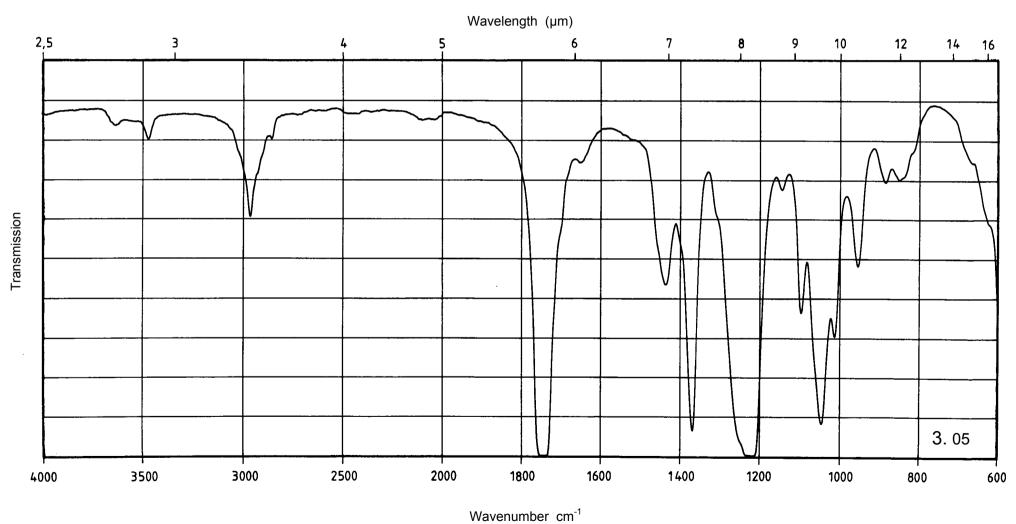


Name	Myristyl lactate	Sample preparation Capillary film, sodium chloride cell
	Ceraphyl [®] 50	Reference Air
		Ordinate 0 - 100 %

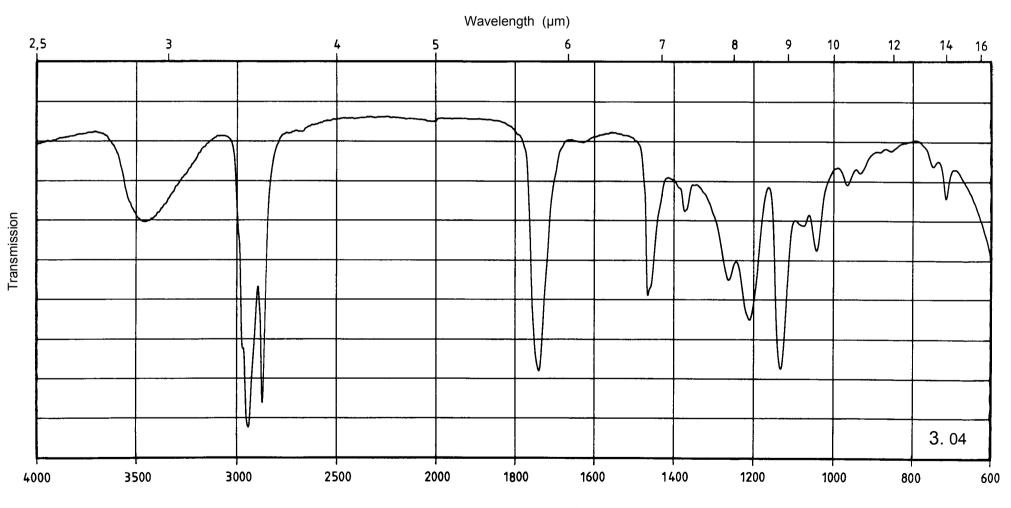


Name	Lauric acid hexyl ester	Sample preparation Capillary film, sodium chloride cell
	Cetiol® A	Reference Air
		Ordinate 0 - 100 %

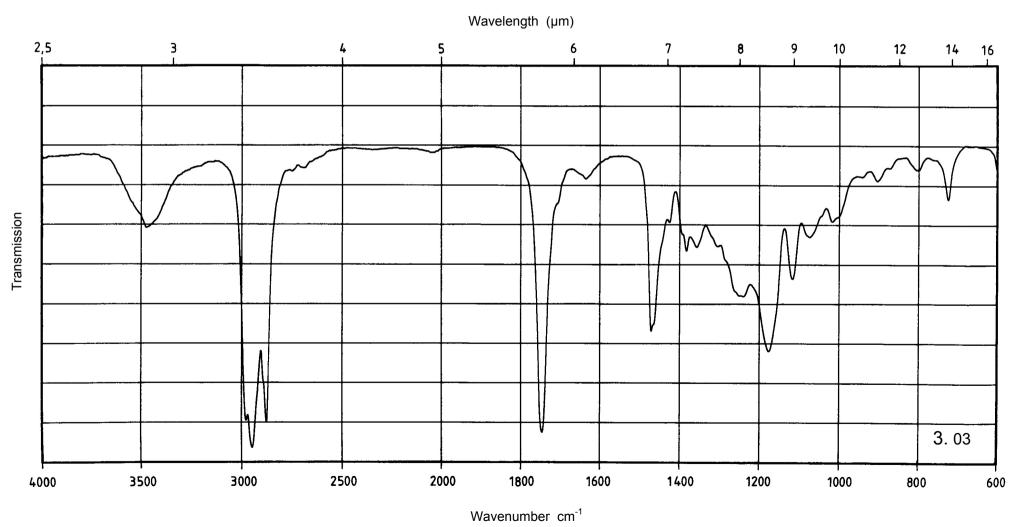


Name	Isopropyl myristate	Sample preparation Capillary film, sodium chloride cell
		Reference
		Air
		Ordinate
		0 - 100 %

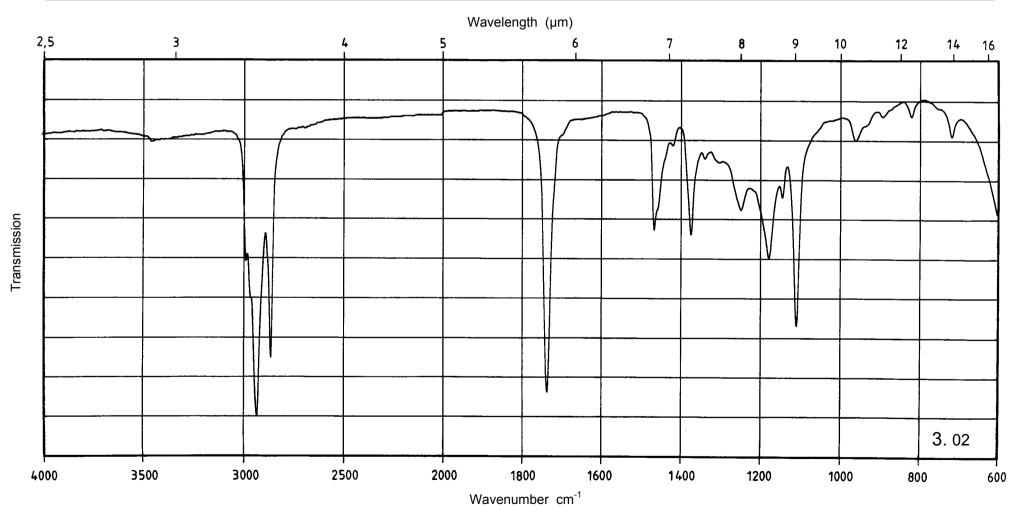


Name	Dehymuls [®] E	Sample preparation Capillary film, sodium chloride cell
		Reference
		Air
		Ordinate
		0 - 100 %

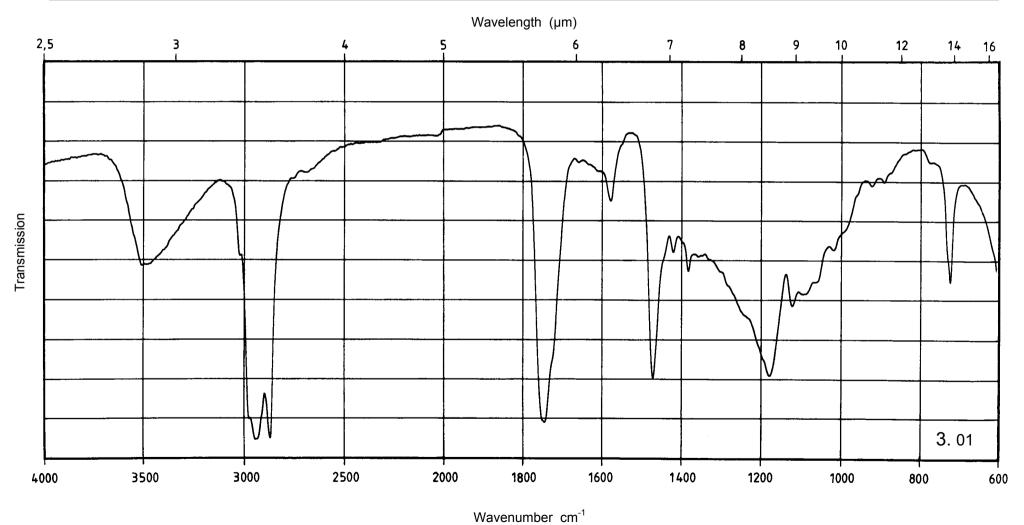


Name	Tris(hydroxymethyl)aminomethane	Sample preparation
		Potassium bromide dispersion
		2 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

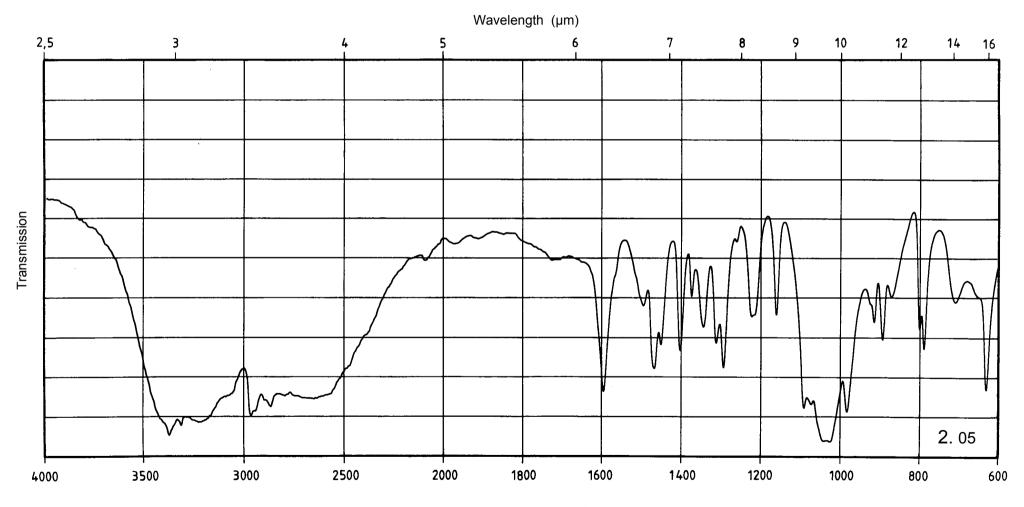


Name Triethanolamine	Sample preparation Capillary film, sodium chloride cell
	Reference Air
	Ordinate 0 - 100 %

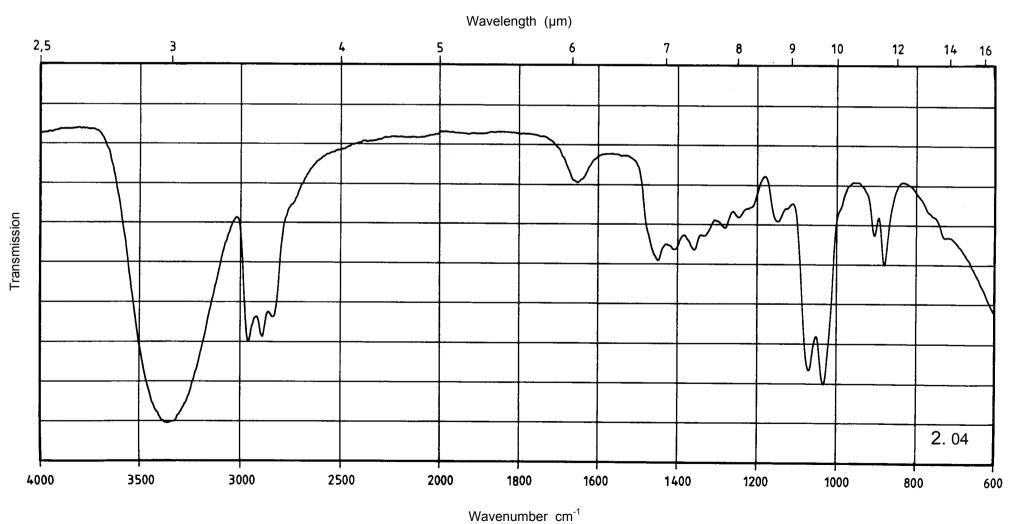


Name	Myristyl picoliniumchloride	Sample preparation
ranic		Potassium bromide dispersion
		1.5 mg / 300 mg
Quartresin [®]	Reference	
	Quartiesin	Air
		Ordinate
		0 - 100 %

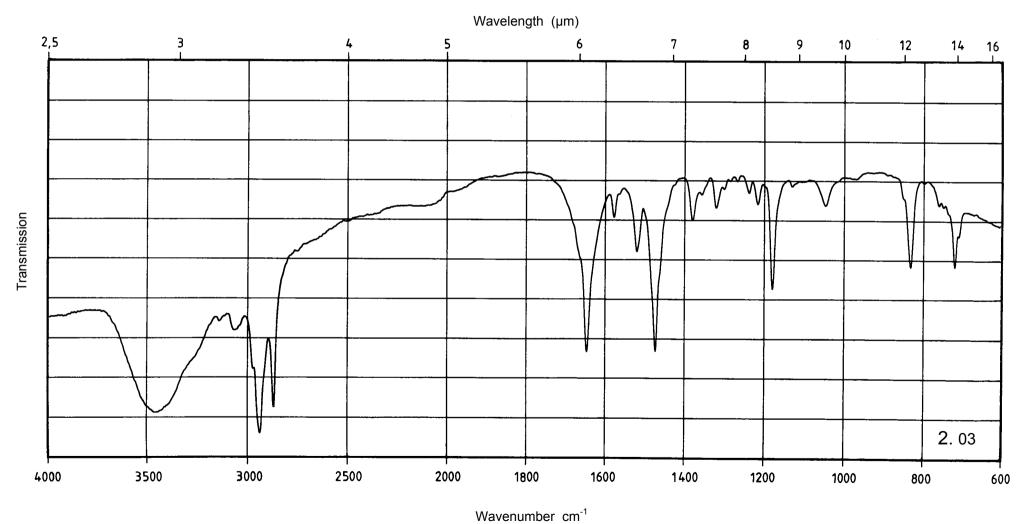


Name	Diethanolamine cetyl phosphate	Sample preparation
Name		Potassium bromide dispersion
		1 mg / 300 mg
	Amphisol [®]	Reference
	Amphisoi	Air
		Ordinate
		0 - 100 %

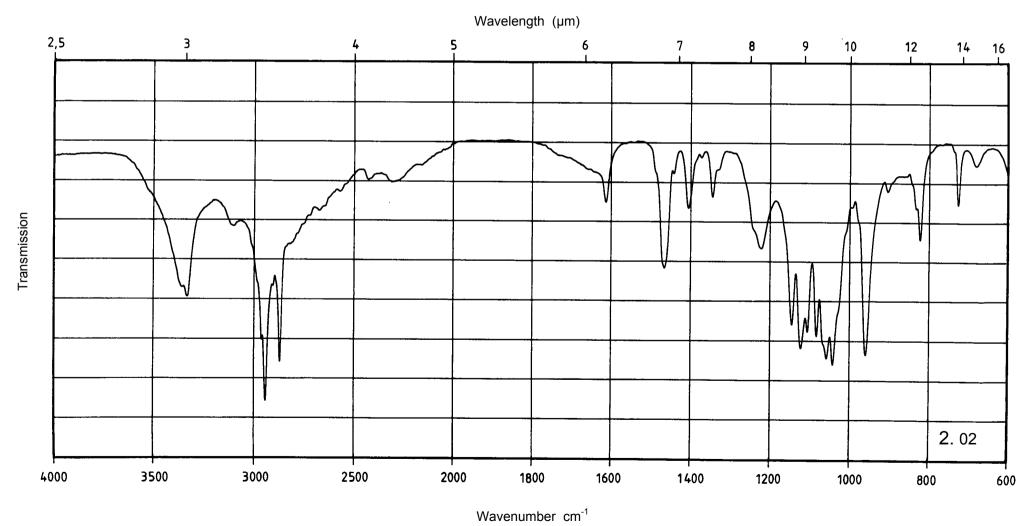


Name	Diethanolamine	Sample preparation Capillary film, sodium chloride cell
		Reference Air
		Ordinate 0 - 100 %

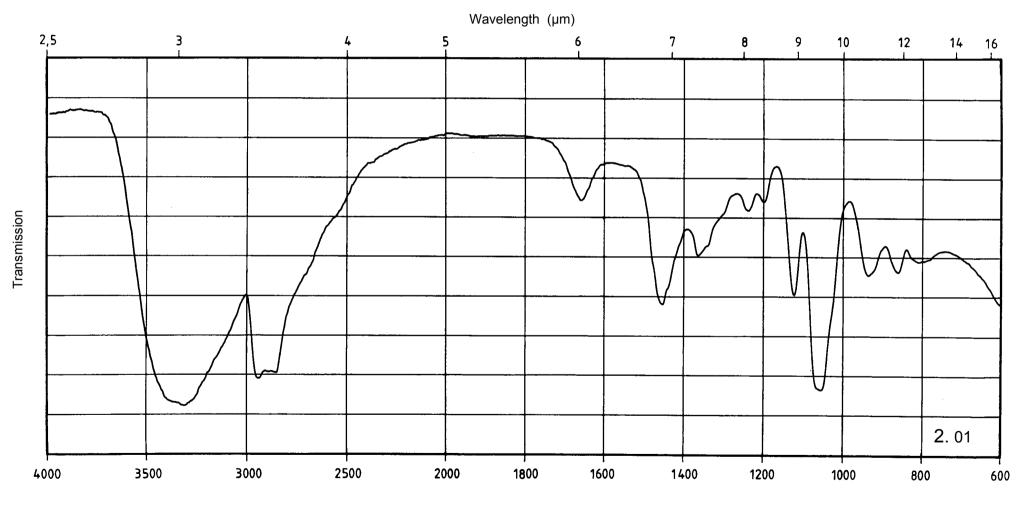


Name	Calcium saccharate	Sample preparation
		Potassium bromide dispersion
		2 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

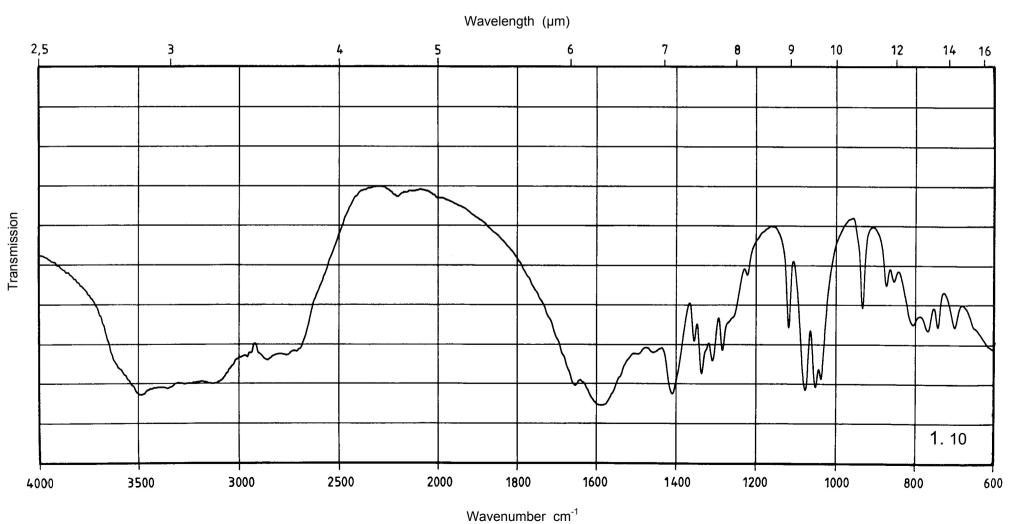


Name	Saccharin sodium, dihydrate	Sample preparation
		Potassium bromide dispersion
		2 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

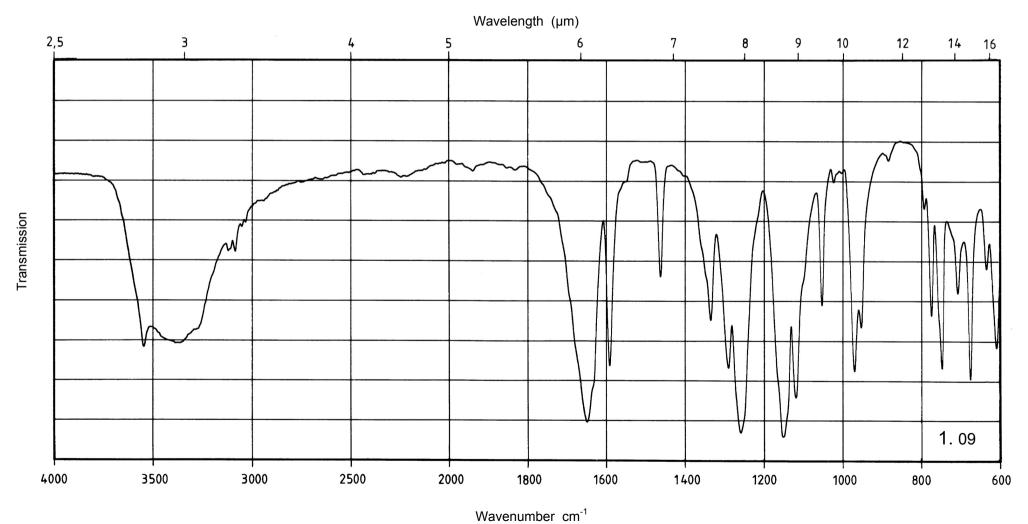


Name	Saccharin	Sample preparation
Ivanic		Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

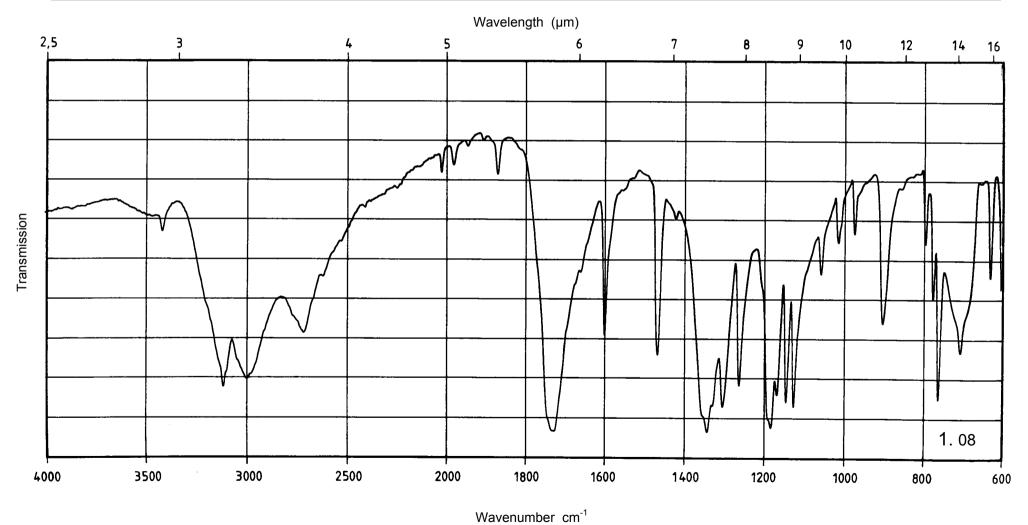


Name	Sodium citrate, dihydrate	Sample preparation
, tamo	Social Trace, and areas	Potassium bromide dispersion
		2 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

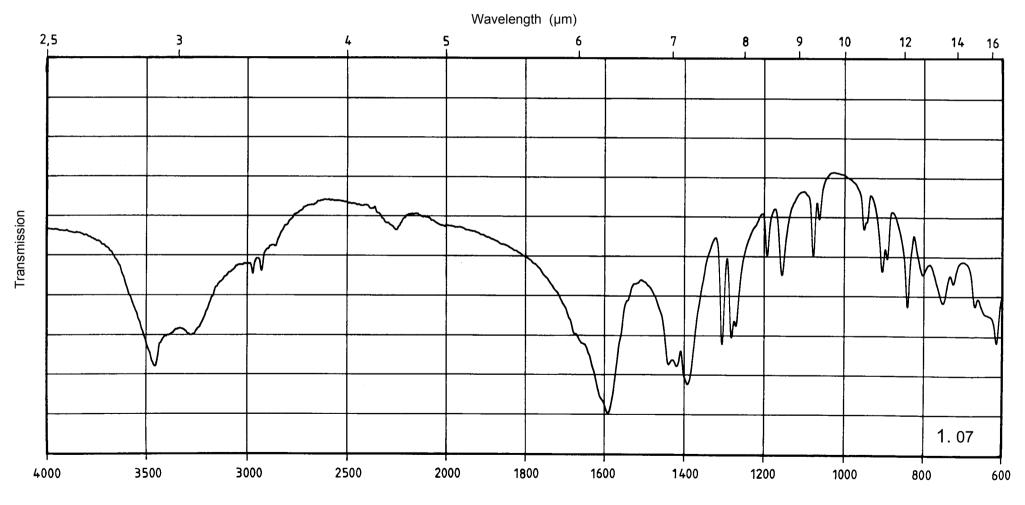


Name	Glycine	Sample preparation
Ivanie	Cryonic	Potassium bromide dispersion
		1.5 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

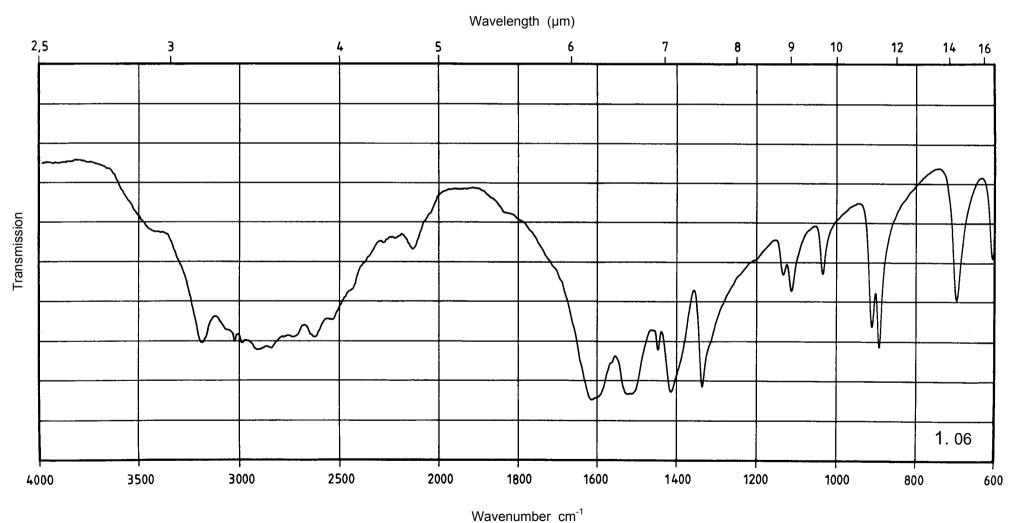


Name	Fumaric acid	Sample preparation		
Ivanic		Potassium bromide dispersion		
		1.5 mg / 300 mg		
		Reference		
		Air		
		Ordinate		
		0 - 100 %		

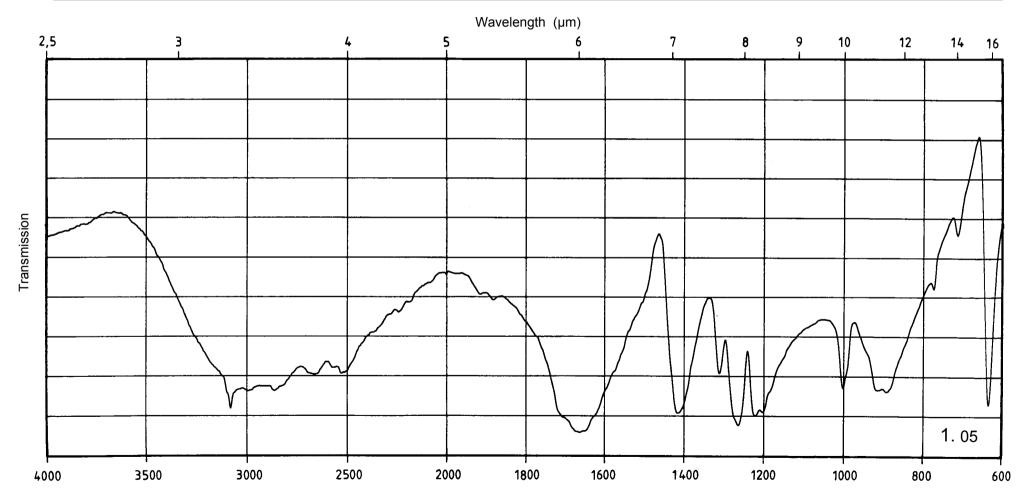


Name	Edetic acid, disodium salt, dihydrate	Sample preparation		
Ivanic	Eactio dola, diocalam sait, amyarate	Potassium bromide dispersion		
		1.5 mg / 300 mg		
	Idranal ® III, Komplexon ® III, Titriplex ® III	Reference		
	idianai ili, Kompiexon ili, himpiex ili	Air		
		Ordinate		
		0 - 100 %		

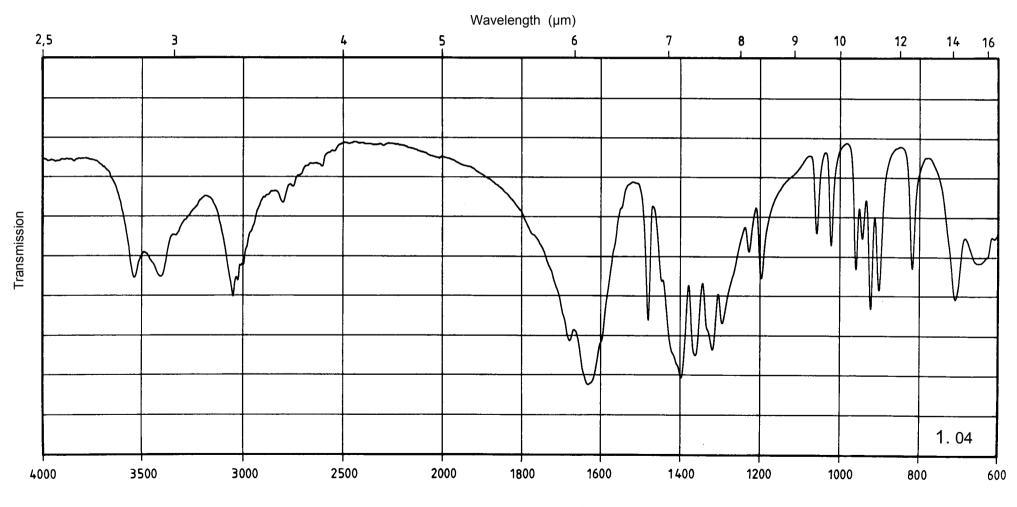


Name	Edetic acid	Sample preparation	
Ivanic	Edetio dold	Potassium bromide dispersion	
		1.5 mg / 300 mg	
	Idranal [®] II, Komplexon [®] II, Titriplex [®] II	Reference	
	iuranai ii, kompiexon ii, nunpiex ii	Air	
		Ordinate	
		0 - 100 %	

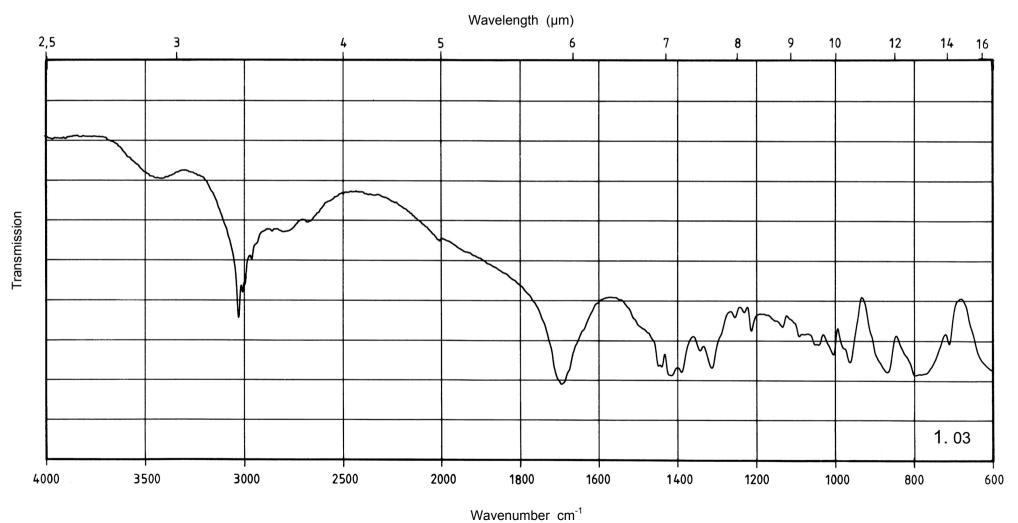

Name	Citric acid, monohydrate	Sample preparation
		Potassium bromide dispersion
		2 mg / 300 mg
		Reference
		Air
		Ordinate
		0 - 100 %

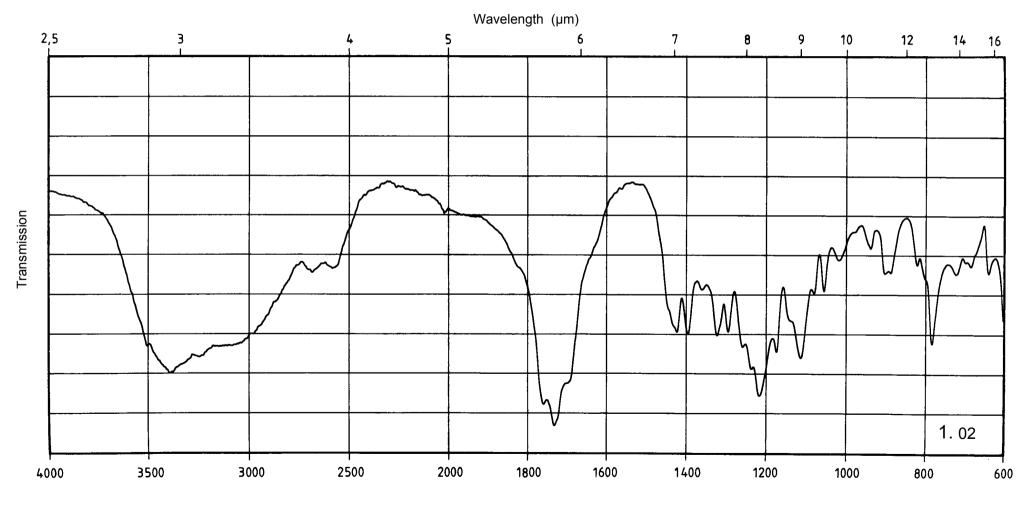

Name	Citric acid, anhydrous	Sample preparation		
Ivanic	Oitho doid, diffrydrodo	Potassium bromide dispersion		
		2 mg / 300 mg		
		Reference		
		Air		
		Ordinate		
		0 - 100 %		

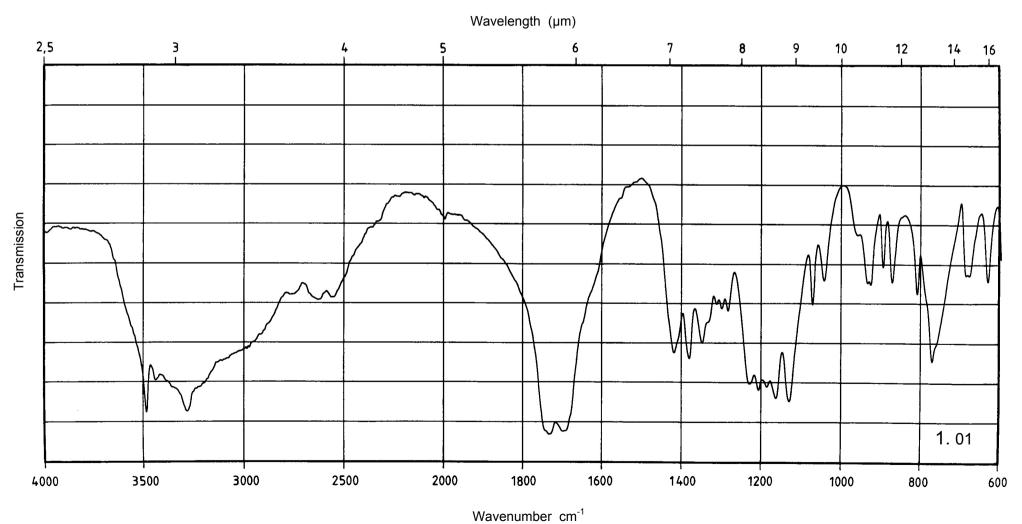

Specfinder UV

increased absorption maximum reduced absorption maximum ↓ ↑↓

increase in one maximum, fall in another

bathochromic shift hypsochromic shift

increase + bathochromic shift fall + hypsochromic shift Χ

considerable change in spectrum

X 283 considerable change in spectrum + bearing a new maximum at 283 nm

 ${\to} \mathsf{Z}$ increasing degradation

Distinct long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1% 1cm	Minima (nm) in methanol	Chang ionis		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
505nm	750	356, <u>235</u>	1315	430, 305			Pyrvinium embonate	217
494nm	219	492, 290, <u>233</u>	675	280, 244		→Z	Daunorubicin-HCI	9108
481nm	207	287, 251	816	368, 270, 216	V		Idarubicin-HCI	9152
480nm	190	337, <u>237</u>	420	400, 297	\downarrow		Rifampicin	420
476nm	215	492, 287, <u>233</u>	668	278, 242	•	→Z	Doxorubicin-HCl	9111
476nm	206	494, 287, <u>233</u>	659	242		<i>→</i> Z	Epirubicin-HCI	9113
462nm	2037	262	2025	220	\downarrow		Acriflavinium Chloride	902
461nm	Z	378, <u>234</u>	2020	422, 320, 224	\rightarrow	Z	Metaclazepam-HCl	2632
440nm	1040	070, <u>201</u>		122, 020, 221		\rightarrow	Ambazone	1220
440nm	380	356, (225)	970	293	←	\uparrow	Carbazochrome dihydrate	3073
440nm	187	240	260	325, 230	- ' -		•	9107
440nm	214		1238		7	کZ لا	Dactinomycin Nitroxoline	1235
440nm 434nm	288	355, <u>238</u>	1160	293, 220		Z	Amsacrine	9130
	223	414, 264, <u>248</u> 269, 224		325, 257, 228 350, 240	Z			
430nm		ļ —	1010		 	<u> </u>	AnhydrotetracyclHCl	2918
430nm	242	373, <u>271</u>	1760	317, 225	· ·	\rightarrow	Ethacridine lactate	913
428nm	491	<u>348</u> , 233	858	288, 217	*	1	Anethole Trithione	3094
426nm	175	<u>271</u> , 224	960	346, 239	<u></u>	<u> </u>	Epianhydrotetracycline-HCl	2921
422nm	420	399, 380, <u>260</u>	3500	410, 386, 224	1	↓	Aminoacridine-HCI	936
408nm	165	<u>292</u> , 228	570	338, 249	←	\rightarrow	Dipyridamole	1410
405nm	22	318, <u>303,</u> 288	660	313, 294, 254			Nystatin	419
395nm	935	(278), <u>238</u>	415	300, 224	\rightarrow	\rightarrow	Phenazopyridine-HCI	123
383nm	574	333, 267, 224		343, 287, 251		X 421	Topotecan-HCI	9156
382nm	291	264	378	315, 242	←	Ľ	Rutoside	2808
381nm	838	309, (226)	395	328, 263, 215	\downarrow	7	Dantrolene sodium salt	2420
380nm	160	<u>296</u>	875	344, 243	\downarrow	Z	Obidoxime Chloride	3058
375nm	635	<u>230</u>		305	Х	\rightarrow	Diminazene aceturate	937
374nm	908	288, 253	610, 446	319, 230		И	Nifurprazine	917
374nm	126	<u>236</u>	702	312, 212			Nilvadipine	1465
373nm	260	<u>268</u> , 232	349	300	↓	Z	Chlortetracycline-HCl	406
373nm	306	<u>268</u> , 228	317	302	Ľ	7	Demeclocycline-HCI	441
369nm	173	300, <u>288</u>	266	312, 296, 269			Pamoic acid	3010
368nm	66	302,290,278, <u>237</u>	1020	313, 211			Benproperine embonate	808
368nm	318	256, (220)	329	296, 234		7	Epitetracycline-HCl	2920
368nm	680	247	393	290	Ľ	7	Sulfasalazine	9109
368nm	327	288, <u>256</u>	376	317, 244			Thiocolchicoside	2433
368nm	775	266	568	315	<u></u>		Triamterene	1313
366nm	903	264	563	306	$\uparrow\downarrow$	\rightarrow	Nitrofural	918
365nm	938	286	535	315, 225	Z	Z	Nifuroxazide	957
364nm	299	<u>266</u>	335	294	1	\rightarrow	Oxytetracycline	410
363nm	139	281, <u>239</u>	524	322, 268		1	Rafoxanide	216
363nm	367	270, (220)	366	298, 235	Ľ	7	Tetracycline	429
363nm	331	270, (220)	345	298, 234	Ľ	7	Tetracycline-HCI	414
362nm	186	<u>237,</u> (221)	524	295, 231			Felodipine	1454
362nm	477			281	Ľ	\downarrow	Floctafenine	124
362nm	801	253, 219		283, 239		7	Olsalazine disodium salt	751
362nm	294	<u>269</u> , 220	349	299		Z	Rolitetracycline	413
361nm	738	<u>286</u>	660	306, 240			Amiloride-HCl	1324
361nm	119	<u>238</u>	342	289, 226			Amlodipine besilate	1460
360nm	418	288, <u>256</u>	467	308, 275		Х	Dithranol	1103
360nm	293	<u>266</u> , 219	364	295, 234	Ľ	7	Oxytetracycline-HCl	430
359nm	1300			265			Etretinate	1110
358nm	522	256, 228	1163	284	\leftarrow	l	Glafenine	107

Distinct long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1% 1cm	Minima (nm) in methanol		ges by ation	Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
357nm	1080	221	672	296	7	7	Iprazochrome	30156
357nm	69	231	778	306	Х	Ľ	Medazepam	2611
357nm	669	269	427	300, 251	1	\rightarrow	Nitrofurantoin	919
356nm	207	276, 235	1350	318, 253, 219	X	-	Proquazone	719
356nm	390	505, 235	1315	430, 305			Pyrvinium embonate	217
355nm	235	(282, 274), <u>237</u>	562	323, (263)	Х	1	Benzbromarone	2505
355nm	840	260	450	297	<u> </u>	X	Furazolidone	944
355nm	657	262	348	300, 218	<u> </u>		Nifuratel	956
355nm	130	238	Z	305, 222	 		Nicardipine-HCI	1464
352nm	174	231	815	303, <210			Nimodipine	1466
352nm	488	<u>257</u> , (223)	950	298, 229	1	\downarrow	Cinoxacin	985
352nm	297	269	368	296, 234	L L	7	Doxycycline-HCI	443
352nm	187	235	742	304		•	Nitrendipine	1456
352nm	274	269	890	300, 240		↑	Riboflavin	2807
352nm	257	255	315	282		X 280	Troxerutin	3018
351nm	988	259	437	302, 226		X 200	Tazarotene	1118
350nm	427	242	775	288			Colchicine	986
350nm	1480	272	773	200			Tretionin	1113
349nm	314	245	451	300, 223		7	Metacycline-HCl	445
348nm	180	237	394	298, 225		7	Meclocycline sulfosalicylate	465
348nm	256	248	385	300, 220	7	7	Minocycline-HCl	447
348nm	90	275	313	308, 256	 		Piretanide	1319
347nm	238	286, (222)	429	315, 258, 215	1	Z	Etofenamate	723
346nm	1280	<u>200</u> , (222)	423	313, 236, 213	<u> </u>			1114
346nm	268	280	367	242 250		7	Isotretionin Mefenamic acid	707
346nm	129	200	307	312, 250 276		//	Viquidil-HCl	017
345nm	305	272	1602	322, 295		Ľ	Aminoquinuride-HCl	903
345nm	128	236	543	280, 218		<u> </u>	Nisoldipine	1449
345nm	81	278, <u>234</u>	1525	302, 255			Quinethazone	1312
344nm	01	330, <u>245</u>	1223	291, 223			Alfuzosin-HCl	2034
344nm	485	273	1072	298, 234	K	L/	Enoxacin	9112
344nm	337	331, 255, 220	828	275	1	<u> </u>	Hydroxychloroguine sulfate	9116
343nm	242	330, <u>246</u>	1235	338, 293, 225	!	V	Bunazosin-HCI	630
343nm	89	(271), <u>236</u>	1530	300			Metolazone	1321
343nm	77	283, <u>253</u> , 216	593	238	*	Ψ	Salmeterol xinafoate	2733
343nm	178	(325), <u>247</u>	1230	294, 224			Terazosin-HCl	627
342nm	965	239	1230	262, 217			Altrenogest	2240
342nm	192	<u>289</u>	987	317, 248	X 255	\downarrow	Niflumic acid	718
342nm	397	270	877	295, 236	7 200	Ψ	Trovafloxacin mesilate	9157
342nm 341nm	183	330, <u>246</u>	1015	339, 294, 227	+		Doxazosin mesilate	631
34 mm	145	235	624	280	1	→	Nifedipine	1418
340nm 340nm	987	235	258	277	 .		Norgestrienone	2207
					\rightarrow	\rightarrow	-	
340nm	277	<u>287</u>	587	317, 255	1	↓	Flufenamic acid	706
340nm	100	270 227	4450	299	7	ı	Tritoqualine	327
339nm	178	279, <u>237</u>	1453	301, 258, 213	Z	\	Azosemide	1333
339nm	942	236	199	260	\rightarrow	\rightarrow	Trenbolone 17β-acetate	1014
338nm	507	329, <u>242</u>	862	272, 226	↓	↓	Dequalinium Chloride	1233
338nm	217	<u>248,</u> 267	949	315		71	Sulfaquinoxaline	922
337nm	159	273, <u>233</u>	1390	295, 248	7	Ľ	Furosemide	1309
336nm	117	<u>270,</u> (224)	552	306, 245	1		Isoxepac	714
335nm	116	(260), <u>222</u>	703	295	<u> </u>	Ľ	Bumetanide	1302
334nm	350	(294), 223	398	265	1	Z	Esculin	3021

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
334nm	6.9	'		281			Thioctic acid	30158
333nm	136	278, <u>234</u>	816	257	И		Quinine-HCI	907
333nm	130	280, <u>234</u>	797	303, 257	<u> </u>		Quinine sulfate	972
333nm	156	280, <u>231</u>	1000	303, 257	<u> </u>		Quinine	971
333nm	271	(234)	1000	270		\rightarrow	Gentisic acid	106
333nm	88	317, 270, 260	219	322, 296, 255		→	Nabumetone	746
332nm	209	285, 242	1760	297, 269, 227			Carazolol	2022
332nm	184	282 282	197	306, 260			Caroverine-HCl	2401
332nm	163	286, <u>242</u>	1415	297, 268, 229			Carvedilol	1462
332nm	163	279, <u>230</u>	1022	244			Quinidine	1402
332nm	142		2230	301	\rightarrow		Moxaverine	
-		<u>240,</u> 319	2230		\rightarrow			2444
332nm	527	200 200 200	0.45	269		\rightarrow	Niclosamide	202
332nm	150	292, 268, <u>226</u>	815	307, 283, 250		\rightarrow	Oxomemazine-HCI	319
331nm		319, <u>280</u>	1060	304, 229		V	Grepafloxacin-HCI	9151
331nm	80	316, <u>271</u> , 262	225	320, 296			Naproxen	121
331nm	233	(342), <u>250</u>	1430	296, 224	↑	א	Prazosin-HCI	611
331nm	252	317, <u>278</u>	955	301, 239			Pefloxacin mesilate	9141
330nm	60	<u>233</u>	786	278	\downarrow	\downarrow	Celiprolol-HCI	2028
330nm	141	249	1644	285, 229	7	Ŋ	Halquinol	1215
330nm	268	254, <u>219</u>	692	279, 244	7	\downarrow	Chloroquine phosphate	908
330nm	180	<u>247</u>	1470	302			Khellin	608
330nm	700			275			Loprazolam mesilate	2630
330nm	143	<u>249</u>	1118	285, 222	\rightarrow		Menadione	2804
330nm	67	271, <u>249</u>	388	292, 258, 230			Phytomenadione	2814
329nm	194	<u>245</u>	2090	276, 224	7	Ŋ	Cloxiquine	973
329nm	496	<u>258</u>	1110	275	←	\rightarrow	Nalidixic acid	941
329nm	170	<u>268,</u> 229	675	290, 243	(↓)	Ľ	4-Chloro-5-sulfamoyl anthranilic acid	2915
328nm	69	<u>235</u>	708	284, 210	\downarrow	\downarrow	Acebutolol-HCl	2002
328nm	111			296	X 247		Cinchocaine-HCl	1714
328nm	1110			260			Mercaptopurine	969
328nm	375	287, (226)	432	307, 270, 218			Sulindac	135
327nm	658			260			Diacetylaminoazotoluene	1102
327nm	75	(278), <u>252</u>	370	235	Ľ	Ľ	Trazodone-HCI	2154
326nm	807	244	1180	275, 213		7	Amrinone	1452
326nm	100	270	449	308, 253	\downarrow		Chlorprothixene	2103
326nm	640	246	435	267, 231	<u> </u>	Z	Cynarine	2511
326nm	74	283, <u>234</u>	2070	321, 309, 248			Dibunafon	826
326nm	Z	200, <u>201</u>	20.0	260		7	Enoximone	1463
326nm	<u> </u>	304, 275, <u>243</u>	925	321, 267, 212			Halofuginone-HBr	974
326nm	286	228	470	270			Isradipine	1469
326nm	578	244	255	278, 220	\downarrow	Ľ	Nizatidine	1625
326nm	65	269, <u>229</u>	714	309, 254	*		Zuclopenthixol-2HCl	2184
325nm	64	270, 230	707	306, 255, 217	+		Clopenthixol-2HCl	2161
325nm	Z	235	101	258	+	Z	Dacarbazine citrate	9134
325nm	280	276	1640	306, 236	+	7	Pipemidic acid	978
325nm	300	276	1795	306, 236	+	7	Pipemidic acid sodium salt	978a
325nm	566	240	368	263, 225	7	7	Piroxicam	727
325nm	143	251	920	276, 225	V		Prothipendyl-HCl	
				· · · · · · · · · · · · · · · · · · ·		↓ ∠	•	314
325nm	527 99	228	490 560	272	Ľ	<u>и</u>	Ranitidine-HCl Bendroflumethiazide	1615 1315
324nm		<u>272</u>		296, 238	1			
324nm	118	272	600	298, 238	1	\	Hydroflumethiazide	1316

Distinct long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
324nm	357	(256), (246)		265		\rightarrow Z	Visnadine	1435
324nm	432	254	400	283, 230			Zomepirac sodium salt dihydr	139
323nm	109	254	1133	293, 234	И	\rightarrow	Clioquinol	910
323nm	31	281, <u>231</u>	2170	319, 248			Sodium Dibunate	825
322nm	867	- , <u></u>		260			Hymecromone	2456
322nm	185			275	\downarrow		Methysergide maleate	3065
322nm	1050			234			Retinol palmitate	2810
321nm	220	298, <u>250</u>	396	308, 276	K		Amoxapine	2159
321nm	440	(214)		256			Carbocromen-HCl	1402
320nm	61	233	1050	300	Z	7	Bromazepam	2613
320nm	87	243	680	290, 230	И	Z	Clotiazepam	2620
320nm	210	229	702	287, < 215		R	Cromoglicic acid disodium salt	2426
320nm	210	298, <u>251</u>	390	307, 277, 242	Х	1	Loxapine	2164
320nm	149	297, <u>249</u>	284	306, 275, 241	Х	1	Loxapine succinate	2165
320nm	111	<u>252</u>	668	277, 226	\rightarrow	\downarrow	Oxypendyl-2HCl	326
320nm	43	<u>268</u>	640	297, 235			Pyrazinamide	999
319nm	12	298, 286, <u>243</u>	262	295, 280, 236	V	\	Fencarbamide napadisilate	2410
319nm	66	289	207	317, 245			Propranolol-HCl	2007
319nm	158			303	↓	\leftarrow	Acemetacin	720
319nm	142	245	580	278, 227		Ŋ	Brimonidine tartrate	30167
319nm	531	244	203	267, 225			Ketorolac tromethamine	753
319nm	64	<u>229</u>	1160	291		7	Lorazepam	2609
319nm	200	289	225	263		Ľ	Meclofenamic acid sodium salt	736
318nm	83	<u>250</u>	266	281, 235	K		Bromhexine-HCI	809
318nm	266	(220)		280		Х	Brotianide	206
318nm	Z	<u>226</u>	Z				Famciclovir	9136
318nm	600	<u>303</u> , 288	660	313, 294, 254			Nystatin	419
318nm	55	<u>265,</u> 235	570	288, 247, 222	↓	↓	Thioproperazine mesilate	2179
318nm	57			305			Proglumetacin dimaleate	738
318nm	100			284	←	7	Xipamide	1314
317nm	316	280	1036	247		Ľ	Ciprofloxacin-HCI	9128
317nm	237	305, <u>262</u>	522	312, 248		Z	Hydralazine sulfate	606
317nm	702	238	813	275		K	Cambendazole	207
317nm	1187	<u>302</u> , 289, 220	1270	310, 295, 250	-		Natamycin	488
316nm	75	<u>255</u>	261	279, 236	7	7	Apalcillin sodium salt	457
316nm	79	<u>271</u> , 227	594	293, 243		<u>↓</u>	Bemetizide	1328
316nm	87	<u>271</u>	625	296, 242		71	Butizide	1317
316nm	97	<u>253</u>	1190	299, 230	7	71	Chiniofon	954
316nm 316nm	80 53	<u>272</u> <u>225</u>	585 824	296, 243	X <u>283</u>	<u> </u>	Cyclopenthiazide Clorazepate, Dipotassium	1318 2601
316nm	 85	228	1300	297		7	Desmethyldiazepam	2916
316nm	306	289, 240	419	302, 264, 229	X <u>283</u>	\rightarrow Z	Flavoxate-HCl	2454
316nm	172	(225)	1 18	299		\rightarrow Z	Indometacin	702
316nm	136	(223) 248	873	274		^	Isothipendyl-HCl	304
316nm	580	228	580	271	<u></u>	Γ	Molsidomine	1417
316nm	55	230	1000	293, 217	''	,,,	Lormetazepam	1536
316nm	237	252	2050	286			Phenothiazine	203
316nm	190	277, 239		299, 255, 222			Rifabutin	4104
316nm	68	<u>263</u>	547	290, 233	Ľ		Thiethylperazine dimalate	349
316nm	73	<u>263</u>	592	292, 238	Ľ		Thiethylperazine dimaleate	350
	67	231	937	300, 215	Z	Z	Camazepam	2606

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
315nm	79	'		300	X <u>282</u>		Diazepam	2604
315nm	208	(283), 241	1490	274, 217	X <u>202</u>	Z	Ethaverine-HCl	1427
315nm	84	(200); <u>211</u>	1100	298	X <u>283</u>	\rightarrow	Oxazepam	2605
315nm	0-1	240		278, 210			Moxaverine-HCI	2445
315nm	77	271, 232	900	298, 248	\rightarrow \downarrow		Periciazine	2129
315nm	81	266	792	293, 248	1	\downarrow	Sulforidazine	2152
-		ļ 			<u> </u>	→		
315nm	129	<u>263</u>	1072	288			Thioridazine	2114
315nm	112	<u>263</u>	936	290, 235			Thioridazine-HCI	2133
315nm	146	<u>251</u>	1900	277, 232	\rightarrow	\rightarrow	Chlorquinaldol	1221
315nm	78	270	592	295, 241	1	7	Cyclothiazide	1307
315nm	980	248	257	263	↓	Х	Hydroxytetracaine- HCI	1706
315nm	99	<u>269</u>	668	293, 240	1	Z	Hydrochlorothiazide	1310
315nm	Z	255	Z	277, 232			Tolmetin	713
314nm	68	270, <u>226</u>	891	293, 245		Z	Polythiazide	1325
314nm	76	<u>230</u>	1080	298, 214	X 284		Temazepam	2619
314nm	200	273, 229		295, 248		\rightarrow Z	Vinpocetine	033
313nm	156			270	И		Ergotamine	2003
313nm	131			270		Z	Ergotamine tartrate	2023
313nm	79	248	282	280, 234	Ľ	Ľ	Ambroxol-HCl	823
313nm	58	<u>229</u>	810	300	Х	Z	Fosazepam	1517
313nm	49	<u>224</u>	700	282		Ľ	Gliquidone	512
313nm	24	<u>280</u> , 270	265	245			Naphazoline nitrate	2715
313nm	51	<u>265,</u> 236	498	288, 247, 222			Pipotiazine palmitate	2178
313nm	190	(240)		268			Lisuride maleate	3064
313nm	61	<u>228</u>	920	300	Х		Prazepam	2612
313nm	145	266, <u>232</u>			Х	\rightarrow Z	Proxymetacaine-HCI	1722
312nm	140			270	Я		Ergocristin	2001
312nm	43	<u>228</u>	742	298	Х		Flurazepam-2HCl	1510
312nm	196			269	\downarrow		Methylergometrine maleate	30108
312nm	110	<u>258</u>	895	280, 226	Ľ	\downarrow	Perphenazine	2148
312nm	78	<u>267</u> , 225	560	288, 241		\downarrow	Trichlormethiazide	1326
311nm	249	241	600	264, 222	1	Ľ	Fenyramidol-HCl	2421
311nm	70	254	266	277, 237			Levobunolol-HCl	2031
311nm	441	<u>277</u> , (213)	487	293, 253	Ľ	Ľ	Metoclopramide-HCI	1611
311nm	392	(230)		262, 218	Ľ	\rightarrow	Ornidazole	958
311nm	244	<u>250</u>	1910	272	\	ΧIJ	Papaverine-HCI	2408
311nm	73	257	620	282, 234	· ·		Prochlorperazine maleate	2131
311nm	94			266	+	Ľ	Noscapine-HCI	805
311nm	109	252	473	279, 224			Physostigmine	1951
311nm	354	230	148	262, 218	← ∠	Z	Tinidazole	967
311nm	77	260	697	282, 225	L L	<u>∠</u>	Trifluperazine-2HCl	
-		<u> </u>			<u> </u>	<u> </u>	· ·	2155
310nm	467	288, (278)	424	296, 252	12		Ethyl Biscoumacetate	3043
310nm	72	<u>260</u>	651	281, 224	L 	Ľ	Fluphenazine-2HCl	2124
310nm	261	298, <u>258</u> , 231	990	304, 290, 235	<u> </u>	1	Halofantrine-HCI	9139
310nm	967	226	242	250	1	↓	Tetracaine-HCI	1709
310nm	506			262	Ľ	\rightarrow	Metronidazole	916
310nm	118	<u>255</u>	1002	277, 224			Chlorpromazine-HCl	2102
310nm	255	<u>245</u>	945	295	1	Х	Chlordiazepoxide-HCl	2602
310nm	509			274	K	7	Mebendazole	209
310nm	178	<u>241</u>	Z	271, 221	7 X	\rightarrow X	Oxyquinoline	1217
310nm	520	232	375	253, 218	Ľ	Z	Oxybuprocaine-HCI	1721

Distinct long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1cm	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
310nm	440	285	407	294, 248		1	Phenprocoumon	3011
310nm	180	(260), 230	830	286, 208			Tiotixene	2153
310nm	365	, , <u> </u>		261	Х	1	Warfarin sodium salt	3013
309nm	Z	241	Z	256, 220	7	\downarrow	Actinoquinol sodium salt	30148
309nm	361	(245)		278	Х	X	Clonazepam	1516
309nm	70	<u>259</u>	621	281, 225	\downarrow		Homofenazine-2HCl	2626
309nm	381	275	464	291, 252	\downarrow	\downarrow	Bromopride	1614
309nm	74	252	266	276, 235			Oxyfedrine-HCl	1431
308nm	945	243	302	250, 235			Benoxaprofen	721
308nm	580	278, <u>230</u>	945	292, 248		7	Ethyl Vanillin	3118
308nm	155	240	345	264, 222			Salicylic acid trichlor- isobutyl ester	728
308nm	190	235	587	263, 221		7	Salsalate	150
308nm	81	<u>256</u>	675	278, 225		\downarrow	Thiopropazate-2HCl	2116
308nm	692	224	692	258	1	Ľ	Dehydracetic acid	3105
308nm	122	<u>245</u>	463	268	1	\rightarrow	Mepyramine maleate	306
308nm	332	<u>252</u>	523	285, 236	Х	\rightarrow X	Flunitrazepam	1511
308nm	91	<u>258</u>	853	280, 224			Triflupromazine-HCI	2156
307nm	167			270			Benzydamine	122
307nm	740	242	295	257	\downarrow	X 325	Hydroxyprocaine-HCI	1712
307nm	79	<u>255</u>	616	278, 221	Ľ	Ľ	Perazin dimalonate	2128
307nm	229	243	710	270, 233		Z	Salacetamide	147
307nm	76	226	1020	294, <210	X		Tetrazepam	2623
307nm	660	<u>277</u>	685	280, 243		Х	Vanillin	3103
306nm	170	239	373	262, 223		\rightarrow	Bornyl Salicylate	30166
306nm	128			271			Bromocriptine mesilate	3039
306nm	148	<u>244,</u> 220	497	274, 230	1		Chloropyramine-HCl	351
306nm	412	<u>246</u>	651	273, 229	<u> </u>	,	Chlorotrianisene	2241
306nm	294	242	494	271, 215	Z	↓	Fluvastatin sodium salt	1813
306nm	380	<u>257</u>	600	290, 240	Х	→366	Nitrazepam	1507
306nm	223	237	480	260, 219		Z	Glycol Salicylate	737
306nm	856	222	906	244	_		Trapidil	30165
306nm	522	(254)	5 40	274	<u> </u>		Triclabendazole	227
306nm	161	245	548	272, 221	7	\rightarrow	Tripelenamine-HCI	315
305nm	132	<u>253</u>	990	276, 222		↓	Mequitazine	330
305nm	112	<u>254</u>	840	276, 220	1		Alimemazine tartrate	317
305nm	72	<u>230</u>	701	289, 221		7	Dosulepin-HCl	2182
305nm	96	227	E06	265	<u> </u>		Labetalol-HCl	2020
305nm 305nm	278	237	586 637	260 240	\	\rightarrow	Methylsalicylate Acenocoumarol	703
	510	282				\rightarrow		3020
305nm 305nm	276 304	317, <u>262</u>	522 339	312, 248 270		Z ↑	Hydralazine sulfate Novobiocin	606 409
		<u>248</u>	ააყ			7		
305nm 304nm	237 Z			274 276			Bithionol	1202 1547
304nm 304nm	234	229	240	260, 222	,	\rightarrow	Lamotrigine Piroctone Olamine	1107
304nm	4.3	<u> </u>	Z 4 U	255	← ∠	\rightarrow X	Methyprylon	1519
304nm	4.3 84	248	217	272, 230	↓ ↓	→	Propafenone-HCl	1433
304nm	215	236	455	260	\downarrow	\rightarrow	Salicylic acid isoamyl ester	705
304nm	100	<u>250</u> <u>253</u>	595	277, 229	Ψ	→	Levomepromazine maleate	2125
304nm	149	264	392	287, 248	\		Methaqualone	1506
304nm	164	265, (218)	404	274, 240	\downarrow		Zeranol	3061
304nm	Z	200, (210)	דטד	242	Ψ	↓ ∠	Zopiclone	1550
50411111				LTL	<u> </u>		Zopicione	1000

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
303nm	177	233	587	264, 214	Ľ	Х	4-Aminophenol	2911
303nm	123	247	863	275, 218	Х		Azaperone	2628
303nm	67	<u>258</u>	634	281, 225	→ →	\downarrow	Chlorambucil	951
303nm	264	200		256	<u></u>	77	Ciclopirox Olamine	909
303nm	50	226	564	277, 219	L L	7	Etiroxate-HCl	1804
303nm	286	220	304	267		7	Fenticlor	9115
303nm	373	267, <u>246</u>	462	278, 256, 230		• • •	Ondansetron-HCI	353
303nm	105	253	860	275, 219	\downarrow		Promazine phosphate	2132
303nm	288	235	532		¥ K	7	· ·	118
303nm	570	262	365	262, 225 233	L	7	Salicylamide Tipprofonia apid	730
				245	1	7	Tiaprofenic acid	
302nm	1200	220	402		↓		Butylaminobenzoic acid	2903
302nm	104	247	343	270, 231	Ľ	L 	Clenbuterol-HCl	2419
302nm	88	246	213	270, 228	1	<u> </u>	Etafenone	1421
302nm	64	261	707	283, 227	<u>↓</u>	<u></u>	Melphalan	976
302nm	473			253	Z	Z	Omeprazole	1627
302nm	500			234		71	Pyrithyldione	2607
302nm	285	<u>234</u>	547	259		Ľ	Salicylic acid	1105
302nm	77	293, 260, <u>225</u>	817	299, 246, 217	4		Sertaconazole nitrate	9154
301nm	209			277		7	Oxyclozanide	213
301nm	233	<u>244</u>	521	270	<u> </u>	7	Buclosamide	906
301nm	268	<u>264</u>	409	287, 246	\	↑↓	Pheniramine 4- aminosalicylate	309
300nm	775	262, <u>238</u>	1760				Carprofen	732
300nm	62.4	274	28.2	279		\downarrow	Glibenclamide	505
300nm	68	<u>235</u>	395	283	1	7	Glymidine sodium salt	509
300nm	246	229	785	272, 220			Oxypertine	2177
300nm	108			246			Proscillaridin	1434
300nm	366	<u>288</u>	370	295, 267			Pyrantel embonate	215
300nm	457	<u>266</u>	628	286, 246	1		Aminosalicylic acid sodium salt	904
300nm	92	<u>250</u>	820	274, 219			Profenamine-HCI	2107
300nm	109	<u>252</u>	922	274, 220			Promethazine-HCl	2108
299nm	567	<u>248</u> , 218	1121	276, 231		Ŋ	Methoxsalen	1116
299nm	107			245			Meproscillarin	1448
299nm	536	(263)		241	Я	X258	Methotrexate	955
299nm	311			256		7	Nimorazole	9119
299nm	515			260	↓	Z	Pyritinol	2151
299nm	210			252	Ľ		Timolol maleate	2017
298nm	62.5	<u>251</u>	439	279, 225	↓		Bamipine-HCI	320
298nm	56	<u>236</u>	611	262, 222		Z	Buspirone-HCI	2631
298nm	265	255, 227		278, 242, <220			Cloprednol	1047
298nm	Z	225		256			Carbimazole	2510
298nm	149			266		7	Hexachlorophene	1231
298nm	59	<u>226</u>	600	275		7	Liothyronine-HCI	2302
298nm	995	227		265	↓	Ľ	Ofloxacin	9120
298nm	1068	241	606	250	1	\rightarrow	Tiabendazole	205
297nm	429			239			Cefotetan	496
297nm	439			248		Ľ	Fentonium Bromide	2440
297nm	63			265			Flecainide acetate	1446
297nm	325			255	71	Я	Ketotifen fumarate	2422
297nm	242	<u>228</u>	410	257	71		Salicylic acid sodium salt	1111
297nm	323	253	633	273, 217	7	7	Acetylaminonitro- propoxybenzene	152
297nm	458	254	278	266, 245			Oxetorone fumarate	3066

length (min) Max. in methanol 297nm 672 223 297nm 160 257, (296nm 437				ation		No.
297nm 672 223 297nm 160 257, (acidic	alkal.		
297nm 160 257, (252	245, 211	\downarrow	Ľ	Procaine phosphate	1720
		276, 245	•		Vincristine sulfate	968
	(==:)	278	\downarrow	7	Albendazole	221
296nm 178 222	1030	252	*	<u></u>	Alizapride-HCI	328
296nm 40 257	349	228	\downarrow		Aprindine-HCl	1436
296nm 65 289, 2		258	\downarrow		Diiodotyrosine	2310
296nm 108	000	274	¥ K		Doxepin-HCl	2110
296nm 32 244	280	276, 237			Ketoconazole	975
	(265, 255) 14.7	276, 251			Normethadone	804a
296nm 570	(203, 233)	268		Γ	Oxfendazole	211
296nm 544 (246)		266	7	7	Oxibendazole	212
296nm 815		243	<u>↓</u>	<u> </u>	Procaine-HCI	1708
296nm 186 263	558	283, 246	Ψ		Zolpidem tartrate	1549
296nm 186 263	556	283, 246			Zotepine	2188
295nm 193 234	325	265, 222	Ľ		Clomifene citrate	2210
295nm 193 259, 2		285, 253, 222	X 240			2190
295nm 356	220	240	X 240	Ľ	Clozapine Penicillin G Procaine	412
	228 1186	247	Ľ.	<u> </u>		30173
				7	Sumatriptan succinate	
295nm 387 268	345	242	\rightarrow		Tienilic acid	1323
294nm 72		267	-	Х	Dembrexine-HCI	3099
294nm 490		270	<u> </u>		Fenbendazole	201
294nm 575		232	\rightarrow	Z	Nitrefazole	3087
294nm 13.3 258	13.5	274	1		Normethadone-HCI	804
294nm 1.6		253			Prasterone	30164
294nm 168 <u>266</u>	274	287, 244			Reserpine	605
294nm 74 (240)		256	7	<u>L</u>	Thalidomide	1531
294nm 125 235	121	274, 255, 226		<u> </u>	Paroxetine-HCI	2187
293nm 373 250	442	270			Brivudine	9131
293nm 146		257		X 311	Bromchlorophene	1229
293nm 50 284	51	(289), 258	1	<u> </u>	Diclofenamide	3077
293nm 23 (283)		235			Flunarizine-2HCl	3045
293nm 39 273, <u>2</u>	<u>253</u> 273	280, 271, 231	\downarrow		Fluspirilene	2143
293nm 740 235	570	261	\leftarrow	1	Uric acid	2906
293nm 13.8 (259)		277			Methadone	110a
293nm 12.3 258	13	274	1		Methadone-HCl	110
293nm 213		276	7		Mirtazapine	2191
293nm 78		272	\downarrow	K	Nomifensine maleate	2106
293nm 466		250	\downarrow	\downarrow	Protriptyline-HCI	2150
293nm 64 257, 2	<u>218</u> 817	273, 237	Ľ	ĸ	Tertatolol-HCl	1468
292nm 66 240	526	274, 220		7	Antazoline-HCl	301
292nm 89 <u>249</u>	262	271, 225	\downarrow	\downarrow	Ajmaline	1423
292nm 1246 220	538	240	↓x	Ľ	Benzocaine	1701
292nm 438 229	361	252	1	Z	Cefixime trihydrate	497
292nm 570 228	520	249	←	\rightarrow	Dipyridamole	1410
292nm 566 (252)		271	*	-	Fendosal	709
292nm 46 283, 2		221		Z	Kawain	2139
292nm 651 (285)		288, 262	<u></u>	7	Parbendazole	214
292nm 236 (320)		247			Propranolol	2007a
292nm 412		255		L	Pyridoxine-HCl	2806
291nm 560 221	595	247	Z		Calcium Folinate Pentahydr.	30152
291nm 139 <u>231</u>	322	249	_	_	Bunitrolol-HCl	2004

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
291nm	200	227	326	253	Ľ	7	Butylmethoxyphenol	3116
291nm	144	232	358	260			Carbamoylphenoxy-acetic acid, sodium salt	749
291nm	56	283	54	264	\	X 258	3,5-Dibromo-4-hydroxy benzenesulfonic acid, sodium salt	1230
291nm	74	282	76	262			Fluorescein dilaurate	30103
291nm	700			225			Medrogestone	2225
291nm	523	254	197	265, 231	↓	1	Methazolamide	1329
291nm	128	227	264	250, 220	\	\downarrow	Midodrine-HCI	2720
291nm	203	<u>281</u> , 271	295	248	\downarrow	\downarrow	Naphazoline-HCl	2715a
291nm	153	280, <u>222</u>	865	289, 244	\downarrow		Pergolide mesilate	30161
291nm	104	225	300	250, 210	Ľ	7	Prenalterol-HCl	2026
291nm	217			265		+	Sildenafil citrate	30172
290nm	549	240	413	254			Befunolol-HCl	1461
290nm	479	223	430	252		Z	Amezinium Metilsulfate	2723
290nm	690	(324)		267			Griseofulvin	407
290nm	210	<u>230</u>	532	258			Ethenzamide	105
290nm	385	285, 245		262	Х		Luxabendazole	229
290nm	129	238	151	259, 222			Piperonyl Butoxide	30162
290nm	76	<u>230</u>	1385	278		Χĸ	Clobazam	2603
290nm	78	<u>226</u>	331	249, 210	\downarrow	\rightarrow	Pramocaine-HCI	1724
290nm	381	248	275	267, 237			Sulbentine	9123
289nm	354	248	267	267, 235	←	1	Bensuldazic acid sodium salt	942
289nm	35			260		7	Buprenorphine-HCI	137
289nm	2,1			227			Camphor	011
289nm	435			254	X	٧	Ethionamide	914
289nm	391			250	Z		Pantoprazole sodium salt	1626
289nm	680			240	\downarrow	¥	Procainamide-HCI	1412
289nm	70			269			Sulpiride	2109
289nm	631			218			Megestrol 17α-acetate	2235
289nm	344			245			Torasemide	1334
288nm	1335			238	$\uparrow\downarrow$	¥	Aminobenzoic acid	2901
288nm	62	279	66	285, 249			Epimestrol	2223
288nm	< 1			256			Erythromycin	423
288nm	181	<u>225</u>	895	249, 215	\downarrow	\downarrow	Metergoline	30160
288nm	42		·	264	Ľ	7	Nalorphine-HBr	3050
288nm	460	278	456	264			Oxantel pamoate	210
288nm	169	<u>265</u>	320	237			Pindolol	2008
288nm	651	259	503	266, 235	\downarrow	Χ	Phthalylsulfathiazole	920
288nm	53	244	168	265, 223			Prajmalium Bitartrate	1453
288nm	835	258	633	227	K	←	Sulfathiazole	933
288nm	202			260	Ľ	↑	Tetroxoprim	966
288nm	67	235	383	267, 228			Tiapride-HCI	3068
288nm	255	233	587	263, 218			Tioxolone	1117
288nm	232			258	K		Trimethoprim	935
288nm	196	<u>243</u>	939	266, 218	←		Piribedil	30163
287nm	145	240	390	260, 227	\downarrow	7	Carbuterol-HCI	2436
287nm	55	280	60	248	Ľ	Ľ	Dextrometorphan-HBr	806
287nm	294	230	270	257		7	Domperidone	338
287nm	465			250	ĸ		Famotidine	1618
287nm	82	224	317	249	\downarrow		Mefexamide-HCI	2166
287nm	61	278	65	275, 238			Mestranol	2213

Distinct	_ 1%	Fronth our many sine o	_ 1%	Minima (nun) in	Chan			Connector
long wave length	E 170	Further maxima (nm) in methanol	E 170	Minima (nm) in methanol	Chang ionis	ation	Pharmaceutical substance	Spectr. No.
Max. in					acidic	alkal.		
methanol	E2			262			Marabina	1110
287nm	53			263 262		7	Morphine LICI	111a
287nm	44	(220)			L L		Morphine-HCI	111
287nm	910	(220)		251	7		Mopidamol	977
287nm	117			256		<i>7</i> 1 →	Pirisudanol dimaleate Silibinin	2174
287nm	446	007	000	253	↓			3012
287nm	792	237	320	252	1	7	Thiopental Sodium	1522
287nm	46,2			257	↓	7	Nalbuphine-HCI	3085
287nm	183			252	1	↓	Nicergoline	025
287nm	731	260	603	269, 229	L	K	Sulfathiazole sodium salt	933a
286nm	423			244	Ľ		Cidofovir	9150
286nm	88			265		L	Crococonazole-HCl	9133
286nm	5.0			246		K	Dehydrocholic acid	1602
286nm	361			255		_	Dimetacrine tartrate	2142
286nm	46			263	<u> </u>	1	Ethylmorphine-HCl	803
286nm	486	(000, 015)		239	1		Fenofibrate	1811
286nm	273	(320, 315)		248			Mebhydrolin napadisilate	325
286nm	61			268		Z	Tolycaine-HCI	1717
286nm	365			260	Ľ	1	Pyrimethamine	979
286nm	52.8	050		264	.,		Codeine	801a
285nm	268	250	200	263, 240	X		Astemizole	343
285nm	39			263	1	↓	Codeine phosphate	801
285nm	338			236			Danazol	2513
285nm	53			259			Dihydrocodeine	802a
285nm	39.8			257			Dihydrocodeine tartrate	802
285nm	94			253			Benzquinamide	333a
285nm	(0.7)						Drostanolone 17β- propionate	2246
285nm	840			230			Dydrogesterone	2247
285nm	39			263			Pholcodine	822
285nm	700	245	296	256, 221		\rightarrow Z	Piprozolin	1621
285nm	640	258	629	268, 224	\downarrow	Х	Succinylsulfathiazole	945
285nm	454	(237)		258			Carbamazepine	1515
285nm	150	<u>249</u>	399	273, 239	Ľ	7	Acetarsol	901
285nm	670			226			Potassium Canrenoate	2508
285nm	244			225		7	Oxymetholone	2227
285nm	741	<u>245</u>	325	257, 222	←	Z	Etozolin	1327
285nm	605	262, <u>231</u>	1650	272, 248		\downarrow	Minoxidil	615
285nm	163	(295), <u>225</u>	1088	246	\	\rightarrow	Naratriptan-HCI	30170
285nm	11.7	277, <u>238</u>	300	219	\downarrow		Oxolamine citrate	816
285nm	657	(260)	533	245		\downarrow	Folic acid	2803
284nm	27	276, 270		243			Benzethonium Chloride	1238
284nm	227			244	И	\downarrow	Idoxuridine	938
284nm	410	223	1345	263			Cyproheptadine-HCl	316
284nm	38			263	Ľ	7	Naltrexone-HCl	30110
284nm	43			253			α-Tocopherol acetate	2812
284nm	129			236	\downarrow	Z	Carbidopa	3040
284nm	186			263	K	7	Pirenzepine-2HCl	1612
284nm	450			220			Cyproterone	2222
284nm	325	(325)		248		Z	Hesperidin	3035
284nm	80			255	Ľ	7	Oxymetazoline-HCl	2721
284nm	89			256			Tetrabenazine mesilate	2137
284nm	353			268	\downarrow	X 246	Febantel	222
284nm	84			256			Benzquinamide-HCl	333

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
284nm	548			216			Chlormadinone 17α-acetate	2229
284nm	45	275, 271, 265, 239	389	280, 223			Citalopram –HBr	2189
284nm	31	<u>275, 232</u>	455	(280), 256, 224			Clofenamide	1330
284nm	71			262		1	Etoposide	990
284nm	Z			243		7	Lansoprazole	1624
284nm	35.6			264	\downarrow	7	Naloxone-HCl	3086
284nm	199	231	254	256, (224)	Ľ	Z	Protokylol-HCl	2726
284nm	37.6			265			Oxycodone	130a
284nm	69.3			262	Z	7	Teniposide	9145
283nm	258			251	\downarrow	7	Benzthiazide	1332
283nm	158	275, <u>251</u>	220	279, 235	Х		Clemizole-HCI	348
283nm	94	275, <u>252</u>	128	280, 238			Clemizole Penicillin	479
283nm	38.6			263	Ľ	Ľ	Oxycodone-HCI	130
283nm	40			265			Hydrocodone	807a
283nm	28			263	\downarrow	\downarrow	Hydrocodone tartrate	807
283nm	38.3			262	↓	\downarrow	Hydrocodone-HCI	820
283nm	592	231	680	261, 222	\downarrow	7	Benzarone	3071
283nm	750			243		Я	Resorantel	228
283nm	55	276	58	246			Bupranolol-HCl	2021
283nm	75	278	75	247	↓	7	Terbutaline sulfate	2431
283nm	165			250	Ľ	\rightarrow Z	Dopamine-HCI	2905
283nm	798			230	\rightarrow	\rightarrow	Canrenone	2504
283nm	132			251	K	\rightarrow	Methyldopa	604
283nm	67	247	133	272, 240			Methanthelinium Bromide	1610
283nm	155	273, (294)	128	243			Naftidrofuryl oxalate	015
283nm	227	<u>227</u>	1260	264, 208			Raubasine	3056
283nm	219	<u>226,</u> (289)	1080	248, <210	Ľ	↓	Yohimbine	3060
283nm	242	274, <u>223</u>	2650	258			Terbinafine-HCl	9155
282nm	413			(215)		_	Cyproterone 17α-acetate	2230
282nm	117	225		250		Z	Isoetarine-HCI	2441
282nm	90	(000)		250	1	_	Levorphanol	129a
282nm	52	(220)		245	1	7	Levorphanol tartrate	129
282nm	118	224	222	251	Ľ	Z	Hexoprenaline sulfate	2427
282nm	61 142	246	118	269, 239	L L	7	Propantheline Bromide	1905
282nm		(070) (000)		254		\rightarrow Z	Corbadrine-HCl	2701
282nm	145	(273), (266)		247	←	V	Pimozide	2130
282nm	121 38	(225)	204	250 242	L L	Х	Isoprenaline sulfate	2713
282nm		276, <u>224</u>	304		Ľ	. 7	Metoprolol tartrate	2011
282nm 282nm	141 162			255 252	L L	\rightarrow Z	Dioxyfedrine Epinephrine	2702
282nm 282nm	95	240	590	265	Z	\rightarrow Z	Acetylsalicylic acid	2703
							anhydride	2917
282nm	129	<u>231</u>	289	254			Emetine-2HCl	912
282nm	425			250	Ľ	7042	Diclofenac sodium salt	711
282nm	326	274	20.0	246	1	7 <u>243</u>	Chlorzoxazone	2402
282nm	27.4	274	29,8	278, 256	L L	X 71	Tizolemide-HCl	1320
282nm	127	228	534	250	<u> </u>		Chlorocresol Dantron	1203 1601
282nm 282nm	435 187	<u>250,</u> (500) 234	890 372	301, 277 265, 222	И	X 7	Benperidol	2141
282nm	11.7	<u>234</u> <u>275</u>	12.3	279, 269	<u> </u>		Dicloxacillin sodium salt	442
282nm	833	213	14.3	239			Fenbufen	717
282nm	109	290	95	246			Gramicidin	431
20211111	108	230	90	240		<u> </u>	Graniiciuiii	401

Distinct long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
282nm	91	(224)		251, 216		\rightarrow Z	Norepinephrine tartrate	2717
282nm	<1	(')		235		, <u> </u>	Mesterelone	2226
282nm	140	(316, 302), 222	1095	246			Mefloquine-HCl	9117
282nm	37	276	35	253			Nabilone	2629
282nm	148	(292), <u>227</u>	1293	250			Rizatriptan benzoate	30171
282nm	17.1	274, 266	30,2	279, 270, 246			Sertraline-HCI	2193
281nm	162	(292), <u>223</u>	879	247	V	V	Cabergoline	2515
281nm	46.6	(202), <u>220</u>	070	248	•		Estradiol 17β-undecylate	2244
281nm	73			249	ĸ	Х	Estriol	2243
281nm	129	225	239	250, 215	_	\rightarrow Z	Isoprenaline-HCI	2713a
281nm	180	223	239	259	\downarrow	7	Triclosan	1236
281nm	142	223	405	253	V V	Z	Dobutamine-HCl	2724
281nm	515	227	405 547	245, 218	↓ ↓	<u>∠</u>		
		221	347	·	V		Indoprofen	733
281nm	180	005	200	249		Z	Norepinephrine	2705
281nm	158	225	386	250		7	Oxatomide	339
281nm	151	000	400	250	א	7	Hexylresorcinol	224
281nm	181	233	432	265, 224		7	Droperidol	2123
281nm	111	(292)		244			Dihydroergocristine	602
281nm	48	(00.4)		256	↓	7	Diamorphine-HCI	127
281nm	100	(291)		244			Dihydroergotamine mesilate	2009
281nm	66	(288)		247	\downarrow	\downarrow	Chlorphenesin carbamate	1101
281nm	69			247		7	Ethinylestradiol	2205
281nm	55	273, <u>240</u>	1177	278, 262,215			Letrozole	2257
281nm	54			245	Ľ	7	Levallorphan tartrate	3002
281nm	129			241	₩	\downarrow	Nebivolol-HCI	2036
281nm	82			247	Ľ	7	Pentazocine	125
281nm	21.1			(259)			Ruscogenin	748
280nm	94			252			Ketobemidone	126a
280nm	87			241	\downarrow	7	Ketobemidone-HCI	126
280nm	240	233	480	264, 218	7		Triprolidine-HCI	331
280nm	239			247	←	Ŋ	Vincamine	800
280nm	97	(219)		250		7	Chloroxylenol	1204
280nm	360			249	7	7	Chlorothiazide	1305
280nm	53			263			Meticillin sodium salt	408
280nm		<u>224</u>	490	214	\downarrow		Tiropramide-HCI	2479
280nm	303	(290)		240			Tryptophan	3006
280nm	311	(289)		240		\rightarrow	Tryptamine-HCI	2909
280nm	420			246	\rightarrow	\rightarrow	Pitofenone-HCI	2409
280nm	75			247		\rightarrow	Estradiol	2202
280nm	53			247		\rightarrow	Estradiol 17β-cipionate	2203
280nm	58			247		•	Estradiol 17β-valerate	2232
280nm	54.5			248			Buserelin acetate	2303
280nm	54 54			247			Gonadorelin	2306
280nm	53.1			245			Goserelin acetate	2311
280nm	64			270			Mianserin-HCI	2140
280nm	915			236	ĽΧ	Ľ	Aminohippuric acid	3009
280nm	103	223, (291)		244	2 /	_	Dihydroergocornin mesilate	619
280nm	94	223, (291)		244			Dihydroergocristine	
							mesilate	602a
280nm	102	223, (291)		244			Dihydroergocryptine mesilate	620
280nm	109	229	262	252, 218		7	Nonivamide	735

Distinct long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1% 1cm	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in					acidic	alkal.		
methanol 280nm	53	(288)		250			Quinestrol	2239
280nm		(200)		237	Ľ	1	Terizidone	
	\rightarrow Z	(074) 040	454	237	<u> </u>	↓		490
280nm	7.2	(271), <u>219</u>	454	054 (000)	1	1	Tioconazole	9110
279nm	170	(289), (273)	100	254, (286)	<u> </u>	↓	Indoramin-HCI	623
279nm	42	226	463	249			Clofibrate	1801
279nm	186			244	1	7	4-Cresol	1219
279nm	78			247		\rightarrow	Estrone	2204
279nm	123	<u>230</u>	335	253			Verapamil-HCI	1416
279nm	66	288	62	285, 249			Epimestrol	2223
279nm	53			249		Z	Metipranolol	2014
279nm	36.6	271	35.5	276, 245			Nadolol	2015
279nm	90	<u>235</u>	261	256, 222			Vetrabutine-HCI	2464
279nm	52	272, (218)	55	275, 241			Toliprolol-HCl	2019
278nm	4.9	(269, 265, 259, 253)					Calcitonin	2512
278nm	67	(228)		256			Gallopamil-HCl	1447
278nm	86	<u>272</u>	95	240			Mephenesin	2429
278nm	68	<u>227</u>	289	246, 213		7	Salbutamol sulfate	2430
278nm	100	224	440	245	\downarrow	7	Fenoterol-HBr	2416
278nm	78	<u>224</u>	397	243	\downarrow	7	Pholedrine sulfate	2710
278nm	335	243	390	256, 230	\downarrow	\downarrow	Butaperazine maleate	2121
278nm	65	287	61	275, 238			Mestranol	2213
278nm	84			247			Butyl hydroxytoluol	3104
278nm	370			242		Z	Orotic acid	2908
278nm	346	345, <u>234</u>	1525	255			Quinethazone	1312
278nm	265			245		1	Dimenhydrinate	323
278nm	521			235	\downarrow	\downarrow	Cefazedone sodium salt	476
277nm	19.7	270	22.3	240			Ketamine-HCI	1535
277nm	85	225	500	246,<210		7	Octopamine-HCI	1441
277nm	71	<u>228</u> , (284)	516	253, 215			Bufexamac	750
277nm	Z	====, (=== 1)	0.0	265	Z	Z	Clavulanic acid potassium salt	3075
277nm	403	237	560	262, 227	7	7	Flucytosine	991
277nm	408	201	300	226	+ ''	71	Formestane	2256
277nm	272			247			Theodrenaline-HCl	029
277nm	206			255	↑	\rightarrow Z		
							Phentolamine-HCI	609
277nm	100	005	222	243	<u> </u>	7	Tyramine-HCI	2709
277nm	51	225	300	244	↓	7	Buphenine-HCI	3023
277nm	13.4	0=1		250		7	Oxytocin	2304
277nm	58	271	64	275, 238	1		Alprenolol-HCl	2012
277nm	155	<u>228</u>	938	261	 	\rightarrow	Phenolphthalein	1603
277nm	75			245	↓	7	Orciprenaline sulfate	2712
277nm	48	(286)			X 240		Clopamide	1306
277nm	7.5	<u>268</u> , 262	9.8	274, 243	1		Clobutinol-HCl	810
277nm	487	<u>252</u>	578	220	↓	Z	4-Hydroxyphenazone	2919
277nm	102	<u>225</u>	565	244		7	Ritodrine-HCI	2465
276nm	390			244	1	7	Guanabenz acetate	622
276nm	30	<u>232</u>	441	258, 214	\downarrow	1	Ciprofibrate	1806
276nm	Z	227	Z	273, 245, 219			Fluoxetine-HCI	2186
276nm	30	269, <u>233</u>	240	274, 259, 230			Ticlopidine-HCl	3069
	62	225	386	244	\downarrow	7	Bamethan sulfate	2714
276nm		1		10.0	†	-	Matter day on a relative a 1101	0700
276nm 276nm	75	225	465	243	\downarrow	7	Methylsynephrine-HCl	2706
	75 65	225 226	465 488	243	\rightarrow	Z	Aspirin	102

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1% 1cm	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
276nm	44	(284)	31				Mefruside	1311
276nm	Z	218	Z	244		7	Meptazinol-HCI	151
276nm	 158			243		7	Thymol	1212
276nm	194			246		7	Resorcinol	1104
276nm	239	227	546	254, 212			Glipizide	511
276nm	83	270, (221)	78	273, 244		Z	Isoxsuprine-HCI	3047
276nm	66	218	209	244	\downarrow	7	Metaraminol tartrate	2719
276nm	957	(215)		246		Ľ	Propylthiouracil	2509
276nm	36	229	236	257, 216			Amoxicillin sodium salt	433a
276nm	400			232		\rightarrow Z	Cefuroxime axetil	417a
276nm	141	234	320	258, 219	Х	, _	Fursultiamine	2821
276nm	24,5	(282)	020	245			Nonoxinol	30112
276nm	65	225	398	246	L L	7	Oxedrine tartrate	2732
276nm	58	(283)		243			Roxatidine acetate hydrochloride	1622
276nm	60			245	+		Leuprorelin acetate	2312
275nm	633			242	V	<u> </u>	Azathioprine	3070
275nm	541			244	1		8-Chlortheophylline	335
275nm	69	(218)		243	<u>'</u>	7	, ,	2725
275nm	250	(210)		245	<u> </u>		Gepefrine tartrate Reproterol-HCl	2450
275nm	232			246		\downarrow	Fenetylline-HCI	2162
						<u> </u>	Pirlindole-HCl	
275nm	303	000	0.4	251	1	71		2171
275nm	20	268	21	273, 250			Estradiol 3,17β- dipropionate	2217
275nm	98			240	↓	7	Etilefrine-HCI	2704
275nm	20	<u>268</u>	20	273, 251			Estradiol 3,17β-diacetate	2216
275nm	106			240	\downarrow	7	Phenylephrine-HCl	2707
275nm	70			249		7	Moxisylyte-HCl	2010
275nm	67	224	412	243	1	7	Synephrine tartrate	2708
275nm	111	(216)		242	\downarrow	7	Norfenefrine-HCI	2711
275nm	517			250	K	Х	Apomorphine-HCI	3001
275nm	29	<u>263</u>	35	243		\rightarrow	Phenoxymethylpenicillin potassium salt	411
275nm	54	(283), <u>227</u>	375	253, 213	\downarrow	\downarrow	Atenolol	2013
275nm	41	268	44	271, 241	\	7	Caroxazone	2138
275nm	12.3	282	11.7	279, 269		1	Dicloxacillin sodium salt	442
275nm	96	225	303	247, 214			Methocarbamol	2423
275nm	20.4	<u>267</u>	24.4	244	↑	Ľ	Nefopam-HCI	133
275nm	46	282, <u>222</u>	323	243			Betaxolol-HCl	2029
275nm	13.7	268	14.0	272, 242			Estramustine 17β-(disodium phosphate)	9114
275nm	80	(281), <u>219</u>	443	245			Gemfibrozil	1812
275nm	318	·		229			Oxabolone cipionate	2251
274nm	146			248	7	7	Buflomedil-HCl	3072
274nm	222			245	1	\downarrow	Cafedrine-HCI	021
274nm		230	510	<220			Estradiol 3-benzoate	2215
274nm	45	<u>224</u>	303	264, 216			Dibutylphthalate	3063
274nm	15	<u>265</u> , 259	21	271, 244	1		Tolpropamine-HCl	311
274nm	410			(240)	7	F	Cefuroxime sodium salt	417
274nm	165			243	Ľ	7	3-Cresol	1216
274nm	212	<u>268</u>	220	(272), 256	1	\downarrow	Mazindol	3014
274nm	372	_ 		232	<u> </u>	-	Articaine-HCI	1705
274nm	122	222	371	243	•		Guaifenesin	2406
∠ (∃000	144		07.1	1-10	<u> </u>	<u> </u>	Cadilonoolii	<u>_</u> +00

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Changionis		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
274nm	54	(283)		280, 270	1	Х	Chlortalidone	1303
274nm	616	(=00)		233	 	L K	Chlorazanil-HCl	1304
274nm	28.4	230	405	260, 209	↑		Cocaine-HCI	1716
274nm	114	254, <u>226</u>	439	(265), 218	'		Pipazetate-HCl	812
274nm	308	254, <u>220</u>	700	246			Propentofylline	034
274nm	54			247			Reboxetine mesilate	2192
274nm	18	232	472	255, 214	\downarrow	X 250	Sotalol-HCI	2016
274nm	81	222	270	245, 213	*	7, 200	Viloxazine-HCl	2157
274nm	32	224	218	260, 216	+		Diethylhexyl phthalate	30101
274nm	56	224	375	262, 216			Diethyl phthalate	30100
273nm	54	267	52	270, 242			Penfluridol	2169
273nm	339	201		233	\rightarrow	И	Azidamfenicol	983
273nm	8.5	266	9.9	200	Z	Z	Cloxacillin sodium salt	989
273nm	372	200	0.0	252	7	_	Cytarabine	988
273nm	186			254	7		Cytidine disodium	
2701111	100			204	''		phosphate	3076
273nm	11.4	(281)		279, 250			Enilconazole	9104
273nm	17.5	(282)		279, 248			Isoconazole	994
273nm	15.5	(281)		279, 247			Isoconazole nitrate	995
273nm	217			246		Z	Etofylline Clofibrate	1808
273nm	344			244			Diprophylline	1406
273nm	342			244			Pentoxifylline	006
273nm	475			244			Caffeine	002
273nm	201			237	Ľ	7	Phenol	1210
273nm	77	223		244			Oxprenolol	2005
273nm	744			229	\downarrow	1	Sulfadoxine	926
273nm	74.3	(280), 227	931	255	↑	\downarrow	Benzoic acid	2902
273nm	190			238	\downarrow	7	2-Cresol	1218
273nm	367			244			Pentifylline	009
273nm	69	280	61	255	\downarrow	\downarrow	Fenoprofen calcium salt	131
273nm	15.4	266	16,0	272, 261	<u> </u>	<u> </u>	Tetryzoline-HCl	2722
273nm	66	279	59	245	1	1	Tramadol-HCl	134
273nm	610	232	1330	249	<u> </u>	·	Acipimox	1810
273nm	39	237	389	215			Ambucetamide-HCl	2452
273nm	40	223	407	253			Bisoprolol fumarate	2030
272nm	8.8	(280)	101	255			Guanfacine-HCl	613
272nm	74	278	62	(276), 238	1		Febuprol	1809
272nm	49	265	40	269, 247	*	*	Maprotiline-HCl	2105
272nm	376	200	10	245			Proxyphylline	2449
272nm	12.1	<u>265,</u> (259)	16.6	270, 245		1	Tranylcypromine sulfate	2134
272nm	297	<u>200</u> , (200)	10.0	235	\rightarrow	\rightarrow	Chloramphenicol	405
272nm	39	<u>266</u>	43	269, 248	7		Lidoflazine	1422
272nm	699	200		234		1	Sulfadimethoxine	924
272nm	297			245	*	7	Methoserpidine	603
272nm	729			227	L L	←	Sulfisomidine	934
272nm	335			233	_	_	Parathion	3019
272nm	46	265	38	268, 245			Benzoctamine-HCI	2617
272nm	Z	200	30	200, 270	\	Z	Cefamandole sodium salt	437
272nm	10.3	<u>280,</u> (264)	10.9	248	*		Miconazole nitrate	948
272nm	185	<u>200</u> , (204)	10.8	243		Х	Pyridostigmine Bromide	1956
272nm	5.8	264, <u>258</u>	9.6	270, 262, 245		^	Quinapril-HCl	633
272nm	212		1150	260, <210			Raubasine-HCl	3057
272nm	495	223	1100	226	Ľ		Sulfaguanole	962
Z1 Z11111	490			220	<u> </u>	\	Juliaguariole	902

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Changionis		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
272nm	154				+	\downarrow	Muzolimine	1331
272nm	200	222		243	1	, L	Yohimbine-HCI	3060a
271nm	439	240	529	215			Ceftriaxone disodium salt	463
271nm	177	230	272	254, 217	\downarrow	Z	Latamoxef disodium salt	455
271nm	284	240	495	222	\downarrow	7	Mofebutazone	726
271nm	532			241		(7)	Theobromine	1414
271nm	401			253	Ľ	Ľ	Acetamino phenazone	2914
271nm	54	(277)		240			Penbutolol sulfate	1409
271nm	18.1	(278)		255			Clonidine-HCl	601
271nm	720			225	\downarrow	Γ.	Sulfafurazole	927
271nm	396			243			Etofylline	1407
271nm	1085			235	\downarrow	\downarrow	Toldimfos sodium salt	3005
271nm	215			234	1	←	Isonicotinic acid	2907
271nm	300			225	←	\rightarrow	Sodium Timerfonate	1209
271nm	739	231	573	247	Х	K	Sulfametoxydiazine	930
271nm	27.4			256	1	Z	Benserazide-HCl	3037
271nm	10.8	(280, 264)		247			Econazole nitrate	946
271nm	20.0	263	21.9	268, 228			Dexfenfluramine-HCl	30154
271nm	20.0	265	21.9	269, 232			Fenfluramine-HCl	3044
271nm	30.2	277	24.4	275, 250	Z		Propicillin potassium salt	449
271nm	800	(215)		244	Ľ	Γ.	Sulfaperin	964
271nm	232			235	Я	Я	Chloramphenicol Succinate sodium salt	478
271nm	38.3	216	115	247			Vindesine sulfate	9159
271nm	920			230	\downarrow	Ľ	Sulfacetamide	921
270nm	442			254	K		Nifenazone	734
270nm	73	(277)		(274), 239			Sulprostone	2307
270nm	434			236		\downarrow	Tegafur	9100
270nm	23.2			256			Trimetazidine-2HCI	1450
270nm	10.0			265			Lofexidine-HCI	614
270nm	45	(276, 260)		240			Phenoxybenzamine- HCI	616
270nm	103			256	1	1	Etacrynic acid	1308
270nm	860			242	~	← X	Sulfamerazine	929
270nm	924			230	\	Ľ	Sulfacarbamide	961
270nm	535			242		7	Aminophylline	1451
270nm	347			241		7	Choline Theophyllinate	1458
270nm	730			242	Z	Χ	Sulfalene	9124
270nm	247			244	┌		Xantinol Nicotinate	30118
269nm	127			264	7		Azatadine dimaleate	332
269nm	43	(276)		239			Phenyltoloxamine citrate	340
269nm	1600	(317)		225	\rightarrow		Rosoxacin	467
269nm	757			242	Х	K	Sulfadimidine sodium salt	925a
269nm	812			241	Х	Ľ	Sulfadimidine	925
269nm	7.6	263, 258		267, 261, 255			Bacampicillin-HCl	450
269nm	840	(248)		223	Ľ	Ľ	Sulfamoxole	939
269nm	501			240		7	Theophylline	1415
269nm	785			236	\downarrow	Ľ	Sulfamethoxazole	940
269nm	332	243	270	251, 232	←	\downarrow	Ramifenazone-HCI	132
269nm	326			241	\downarrow		Metixene-HCI	3049
269nm	8.8	<u>259</u> , 254	11.2	229			Phenprobamate	2424
269nm	770			233	И	Ľ	Sulfametrole	963
269nm	652			226	Х	K	Sulfaphenazole	965
269nm	610			240	\downarrow	Γ.	Bentiromide	3095

Distinct long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1% 1cm	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
269nm	34	275	32	255			Carindacillin sodium salt	493
268nm	229	222	990	249	V	V	Bopindolol malonate	2035
268nm	851			240	Ľ	Х	Sulfadiazine	923
268nm	11.1	262, 256	15.1	266, 259, 228		\downarrow	Pargyline-HCl	2111
268nm	220	274	212	(272), 256	1	+	Mazindol	3014
268nm	356	234	408	243	K	←	Dimethylamino phenazone	103
268nm	11.1	(262)		(266)		`	Pemoline	010
268nm	5.6	262, <u>257</u> , 251	8.9	266, 260, 244			Phenamacide	2483
268nm	775		0.0	240	Х	←	Sulfamethoxypyridazine	931
268nm	740			227	1	Ľ	Carbutamide	503
268nm	44	275	41	271, 241	 	7	Caroxazone	2138
268nm	9.8	277, (262)	7.5	274, 243	*		Clobutinol-HCl	810
268nm	373	234	366	237, 215	\downarrow		Clofezone	741
268nm	22	(274)	000	262	T		Flucloxacillin sodium salt	483
268nm	3.1	264, 258, 253		266, 263, 256			Trandolapril	634
267nm	328	(287)		240	+		Mepindolol sulfate	2024
267nm	664	(201)		230	V	V	Urapidil	618
267nm	4.6	261, 256		265, 259, 254		Z	Ampicillin trihydrate	401
267nm	294	238	489	255, 221	1		Phenazone comp. with	
20711111	234	230	409	255, 221	'		salicylic acid	143
267nm	19	<u>236</u>	510	225			Clofenotane	3022
267nm	21	229	309	253	1		Meclozine-2HCl	305
267nm	12.4	275, 260		245		1	Baclofen	2414
267nm	9.2	(274)		254			Butanilicaine citrate	1702
267nm	11.3	(274)		254			Butanilicaine-HCI	1703
267nm	9.2	(274)		254			Butanilicaine phosphate	1704
267nm	4.35	<u>260</u> , 254	5.70	265, 257, 252			Bornaprine-HCI	3038
267nm	492	243	504	255, 233	X 228		Phenazone	115
267nm	180	<u>237</u>	315	260, 216	\downarrow	Z	Cefoxitin sodium salt	439
267nm	18	(264, 257)					Enalapril maleate	626
267nm	248	<u>233</u>	404	251, 220	Х		Thiamine nitrate	2809a
266nm	21.0	<u>261</u>	23.4	264, 228			Fluconazole	9138
266nm	6.8	260	9.1	264, 257			Isopropamide Iodide	2428
266nm	398			229	И	Ŋ	Nifenalol-HCl	2025
266nm	9.75			258	1	Х	Xylometazoline-HCl	2716
266nm	297			235		Z	Dilazep-2HCl	1419
266nm	374			234		\downarrow	Edoxudine	9135
266nm	303			235		Z	Hexobendine-2HCl	1420
266nm	6.7	259, (274)	5.8	242			Clofenciclan-HCl	2115
266nm	551			233	Ľ		Hexamidine isetionate	1237
266nm	47	(273)		261		Χ	Indanazoline-HCI	2730
266nm	509			231	Ľ	Ľ	Pentamidine isetionate	9142
265nm	12.3			242	↓	V	Bambuterol-HCl	2467
265nm	4.5	<u>258</u> , 252	6.2	262, 255, 249			Butetamate-HCI	2470
265nm	7.9	<u>258</u> , 253	10.5	263, 255, 244			Butinoline phosphate	2471
265nm	537			232		Ŋ	Fluorouracil	992
265nm	22.5	(272)		269, 251	1		Thiamphenicol	469
265nm	25.3	(273)		271, 262			Tramazoline-HCl	2728
265nm	Z			239		Z	Cefsulodin sodium salt	453
265nm	11	259	13	263, 257			Lorcainide-HCI	1440
265nm	224	234	265	256	K		Metamizole magnesium salt	138
265nm	6.4	<u>259</u> , 253	8.2	264, 257, 252			Pethidine	113a
265nm	71			248	X		Saccharin	3102

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
265nm	10.4	227, 259	637	248	X	Х	Chlormezanone	2413
265nm	20.4	276, <u>231</u>	575	212	,	<u></u>	Chlorpropamide	504
265nm	4.38	259, 253	5.62	263, 256, 251		•	Adiphenine-HCI	2418
265nm	4.6	259	7.0	263, 257		Z	Ampicillin sodium salt	434
265nm	7.4	258	10.8	263, 256			Azlocillin sodium salt	435
265nm	21.9	271	20.0	269, 232			Fenfluramine-HCI	3044
265nm	12.7	259	17.1		↓		Fedrilate maleate	814
265nm	169			232	7	Z	Iproniazid phosphate	2144
265nm	3.17	259	4.72	263, 257	Z	Z	Mezlocillin sodium salt	446
265nm	240	234	288	254	K		Metamizole Sodium	109
265nm	683			230	\downarrow	Ľ	Sulfacetamide sodium salt	921a
265nm	387			235	•		Zidovudine	9149
265nm	4.6	258	6.6				Benzylpenicillin potass. salt	402
264nm	13.2	<u>258,</u> 252	16.5	261, 254, 248			Benzyl mandelate	2468
264nm	3.1	<u>258</u> , 253	4.1	262, 255, 249			Butetamate citrate	2469
264nm	91	222	500	246, <210	\downarrow		Etofibrate	1807
264nm	15.5	(273)		256			Tocainide-HCI	1444
264nm	440	(2.3)		228			Calcifediol	2818
264nm	206	231	390	(252, 216)	\	Х	Cefadroxil monohydrate	460
264nm	12.0	(258)		243	*		Orphenadrine citrate	2446
264nm	7.2	258, 252		263, 256, 250			Oxetacaine-HCI	1723
264nm	11.4	(272)		248, 269			Praziquantel	226
264nm	8.6	<u>258,</u> 251	11.0	262, 255, 249			Primidone	1521
264nm	879	200, 201	11.0	230	+	Ľ	Sulfaguanidine	928
264nm	265	230	350	248, 220	X		Metyrapone	1322
264nm	5.9	258, 252	7.4	262, 254, 249	\ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \		Methylphenidate-HCl	2113
264nm	14.5	272, 258	7.7	270, 248		↑	Ibuprofen	708
264nm	4.6	257	6.8	263, 256		ı	Benzylpenicillin sodium salt	428
264nm	19.7	228, 275	434	253, 212			Gliclazide	510
264nm	11.1	220, 213	707	236			Mexiletine-HCl	1430
264nm	7.5	258	11.4	262, 257			Tolazoline-HCl	2018
264nm	430	200	11	227			Alfacalcidol	2819
264nm	773			225	7	Ľ	Deanol 4-acetamido-	
20						_	benzoate	2181
264nm	13.8	<u>258</u> , 252	16.5	261, 254, 249			Tropalpin-HCl	2480
263nm	17.1	270, 258		253			Bietamiverine-2HCl	2434
263nm	417	<u>230</u> , (324)	855	257, <210	1	Х	Amiphenazole-HCI	019
263nm	14.8	(271)					Lidocaine	1707a
263nm	14.2	(271)		(257)			Lidocaine-HCI	1707
263nm	13,7	(272)		261			Etidocaine-HCI	1718
263nm	20.8	274, <u>228</u>	453	273, 213			Tolazamide	514
263nm	14.2	(271)		257	↑	↑	Mepivacaine	1710
263nm	184			249	1		Nikethamide	004
263nm	16.6	274	12.2	272			Glibornuride	506
263nm	22.2	274, <u>228</u>	496	271			Tolbutamide	508
263nm	225	(270)		239	1		Nicotinic acid	1802
263nm	479	210		225			Colecalciferol	2811
263nm	11.5	(272)		(258)			Bupivacaine-HCI	1715
263nm	190	·		244	Ľ	Z	Cefradin monohydr./sod. carb.	440
263nm	220	<u>236</u>	371	(257)			Cefalotin sodium salt	415
263nm		<u>235</u>		248, 219	\downarrow	1	Aztreonam	474
	13.2			235			Debrisoquine sulfate	625
263nm	13.2	(271)		233			Debrisoquine sullate	025

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
263nm	201	220		242	1		Nicofuranose	30111
262nm	131			247	<u> </u>		Brompheniramine maleate	329
262nm	259	219	290	241			Suloctidil	3067
262nm	485	2.0		235	←	←	Zimelidine-2HCl	2172
262nm	286	238	352	256, 213	<u> </u>	X	Thiamine-HCl	2809
262nm	152	(256)		240	†		Nicotinic acid benzyl ester	704
262nm	226	(200)	239	2.0			Cefalexin	403
262nm	15.1	268, 256		266, 259, 228		\downarrow	Pargyline-HCl	2111
262nm	182	(269)		235	↑		Pheniramine	312
262nm	1116	(/		225	<u> </u>	Ľ	Sulfanilamide	932
262nm	334			229	↑	\rightarrow	Isoniazid	915
262nm	221			245	<u> </u>	,	Nicotinamide	2805
262nm	138			240	<u> </u>		Pheniramine maleate	308
262nm	433			220	<u>'</u>	7	Acetazolamide	1301
262nm	140			244	1		Chlorphenamine maleate	303
262nm	Z			232	\rightarrow	Z	Cefapirin sodium salt	438
262nm	309			238		Z	Loracarbef monohydrate	4100
262nm	256			228	↑	_	Pyridylcarbinol	3055
262nm	212			250	<u> </u>		Hydroxymethyl-	
20211111	212			200	'		nicotinamide	30122
262nm	286	<u>221</u>	695	238			Mebeverine-HCI	2458
262nm	134			240	1		Nicoboxil	745
261nm	182			236	1		Disopyramide phosphate	1438
261nm	264			226	1		Nicotinyl Alcohol	027
261nm	113			230	1		Nicotinyl Alcohol Tartrate	028
261nm	136			245	1		Carbinoxamine maleate	302
261nm	157			246	1		Bisacodyl	1604
261nm	181			231	1		Nicotine	3003
261nm	251			247			Sulfur	1106
261nm	109			236	1		Doxylamine succinate	2618
261nm	732			230	V	Z	Leflunomide	752
261nm	18.5	(267)		240	\downarrow	Z	Neostigmine Bromide	1955
261nm	23	<u>224</u>	725	244			Tosylamine	3008
261nm	340			258		\downarrow	Trifluridine	9127
261nm	13.1			236	1		Tulobuterol-HCl	2463
260nm	15	254, <u>225</u>	341	256, 244, 217		1	Cyclizine-HCl	352
260nm	680	<u>251</u> , 243	1019	257, 246, <210			Dihydrotachysterol	2822
260nm	48.5			234	1		Indinavir sulfate ethanolate	9153
260nm	13.1	(269, 265, 255)					Metamfetamine	2126a
260nm	14.6			245			Pramiverine-HCI	2448
260nm	12.8			(257)	1		Hydroxyzine-2HCl	2610
260nm	233			239			Phenindamine tartrate	307
260nm	22			250	1		Clotrimazole	911
260nm	5.70	267, 254		265, 257, 252			Bornaprine-HCI	3038
260nm	126			234			Cetylpyridinium Chloride	1225
260nm	7.0			(257)			Dextromoramide tartrate	108
260nm	449			233	K	7	Etifoxine	2146
260nm	12.8	275, 266, 254		246			Loperamide-HCI	1609
260nm	163	(266)		225	1		Betahistine-2HCl	344
260nm	104	(266)		225	1		Betahistine dimesilate	345
259nm	240			232		Z	Cefotiam-2HCI	452
259nm	8.6	265, 253		263, 256, 250			Isoaminile	815

Distinct	40/		40/					
long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1% 1cm	Minima (nm) in methanol	Chang ionis		Pharmaceutical substance	Spectr. No.
Max. in					acidic	alkal.		
methanol 259nm	14.5	(254)		246			Pipradol-HCl	2173
259nm	10.3	265, 253		264, 255, 250			Trospium Chloride	2432
259nm		265, 253					Butamirate citrate	817
259nm	3.83	· · · · · · · · · · · · · · · · · · ·		263, 256, 250 239			Clofedanol-HCl	
259nm	13.0	(265)		245				819 337
259nm	112	264, 253		246			Diphenylpyraline-HCl Propiram fumarate	
259nm	112	(264, 264)		249			Quinaprilat	146 636
259nm	1360	(204, 204)		217		L L	Thiamazole	30174
259nm	14.3			245		Ľ	Terfenadine	342
259nm	5.26	265, 253		263, 256, 250			Valethamate Bromide	2451
259nm	14.0	(268, 252)		239			Fenpiprane-HCI	2417
259nm	11.3	(264, 253)		(250)			Fenpipramide-HCl monohydr.	2417
259nm	566	(204, 255)		227			Adenosine	1401
259nm 259nm	414			232			Adenosine phosphate	30149
259nm	247			231			Adenosine triphosphate	30150
259nm	5.62	265, 253		263, 256, 251			Adiphenine-HCI	2418
259nm	615	200, 200		222	L L	Ľ	Chlorhexidine diacetate	1226
259nm	16.0	(275, 264, 254)		249	↑	L	Chlorphenoxamine	321
259nm	13.6	(275, 264, 254)		249	<u> </u>		Chlorphenoxamine-HCI	321
		L			ļ ļ		·	
259nm	10.6	265, 253		263, 255, 251			Clidinium Bromide	1954
259nm	3.49	265, 253		262, 256, 250			Oxeladin citrate	811
259nm	11.2	269, 254	04	229			Phenprobamate	2424
259nm	723	(320)	21	245			Tolnaftate	947
259nm	13.0	264, 253	000	263, 256, 231			Tribenoside	3059
259nm	574		228				Vidarabine	9147
258nm	1556	(204 252)	226	220			Anethole	824
258nm	10.9	(264, 252)		228			Ephedrine Alfentanil-HCl	003a
258nm	5.8	202 252		044				148
258nm	11.0	263, 253		244			Benzatropine mesilate	1908
258nm 258nm	8.9 4.93	264 264, 252		(263, 255), 238	Z		Benzathine Benzylpenicillin Ciclonium Bromide	495 2437
258nm	7.72	264, 252		262, 255, 250			Cyclandelate	022
258nm	12.8	264, 253		244			Difenidol-HCl	336
258nm	12.9	264, 252		(255, 261), 246			Diphenoxylate-HCl	1617
258nm	5.77	264, 252		262, 255, 249			Drofenine-HCl	2439
258nm	6.72	264, 252		262, 255, 249			Methenamine mandelate	1224a
258nm	3.78	264, 252		262, 255, 249			Pentoxyverine citrate	821
258nm	10.2	264, 253		245			Piritramide	149
258nm	715	204, 200		222	+	Ľ	Proguanil-HCI	9143
258nm	7.2	264, 252		228	*	_	Programme-HCI	026
258nm	8.7	(264, 252)		228			Methylephedrine-HCl	026
258nm	15.0	(263, 253)		238			Prenylamine lactate	1411
258nm	9.6	(264, 252)		228			Metamfetamine-HCl	2126
258nm	11.8	(263, 252)		238			Phenylacetate sodium salt	3004
258nm	8.06	(263, 252)		(255, 250)			Fentanyl citrate	128
258nm	17.1	(263, 252)		228			Benzyl Alcohol	1201
258nm	6.15	(268, 262, 252)		247	1		Camylofine-2HCl	2415
258nm	15.5	(263, 252)		243	'		Diphenhydramine-HCl	310
258nm	9.7	(264, 252)		228			Amfetamine sulfate	2120
258nm	14.7	<u>227</u> , (263)	420	254	↑		Bromazine-HCl	313
258nm	15.9	(252)	1 20	240	'		Polycaine-HCl	1711
258nm	7.4	(264, 252)		262, 254, 249			Methylphenidate-HCl	2113
258nm	12.9	(263, 252)		261, 255, 248			Mandelic acid	3017
2301111	14.9	(200, 202)		201, 200, 240			Ivialiuelle aciu	3017

Max. in methanol acidic alkal. 258nm 16.2 (263, 252) 261, 255, 248 Mandelic acid be 258nm 5.56 (264, 252) 234 Biperiden-HCI 258nm 5.6 (264, 252) 262, 255, 248 Homatropine-HBi 258nm 695 240 ∠ ∠ Chlorhexidine-2H 258nm 370 236 ↑ ↑ Dimetindene mal 258nm 7.25 (264, 253) 263, 255, 228 Mefenorex-HCI 258nm 18.4 (264, 252, 248) (257, 244) Mesuximide 258nm 9.82 (264, 252) 262, 254, 229 Norpseudoephed 258nm 5.51 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 10.6 (264, 252) 262, 255, 250 Phenglutarimide- 258nm 5.78 264, 252 262, 255, 250 Phenglutarimide- 258nm 16.8 238 Fendiline-HCI 258nm 5.64 (264, 252) 230 Fenproporex-HCi <th>drine-HCl -HCl -HCl -HCl -HCl -HCl -HCl -HCl</th> <th>2407 2411 1902 1227 324 3048 1530 3051 3052 3053 1910 032 1455 30104 2455</th>	drine-HCl -HCl -HCl -HCl -HCl -HCl -HCl -HCl	2407 2411 1902 1227 324 3048 1530 3051 3052 3053 1910 032 1455 30104 2455
258nm 16.2 (263, 252) 261, 255, 248 Mandelic acid be 258nm 5.56 (264, 252) 234 Biperiden-HCl 258nm 5.6 (264, 252) 262, 255, 248 Homatropine-HBi 258nm 695 240 ८ ८ Chlorhexidine-2H 258nm 370 236 ↑ Dimetindene mal 258nm 7.25 (264, 253) 263, 255, 228 Mefenorex-HCl 258nm 18.4 (264, 252, 248) (257, 244) Mesuximide 258nm 9.82 (264, 252) 262, 254, 229 Norpseudoephed 258nm 5.51 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 10.6 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 5.78 264, 252 262, 255, 250 Phenglutarimide-Cyclopentolate-H 258nm 9.91 (264, 253) 249 Doxapram-HCl 258nm 8.0 (263, 252) 230 Fengline-HCl 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264,	drine-HCl -HCl -HCl -HCl -HCl -HCl -HCl -HCl	2411 1902 1227 324 3048 1530 3051 3052 3053 1910 032 1455 30104
258nm 5.6 (264, 252) 262, 255, 248 Homatropine-HBI 258nm 695 240 ∠ Chlorhexidine-2H 258nm 370 236 ↑ Dimetindene mal 258nm 7.25 (264, 253) 263, 255, 228 Mefenorex-HCI 258nm 18.4 (264, 252, 248) (257, 244) Mesuximide 258nm 9.82 (264, 252) 262, 254, 229 Norpseudoephed 258nm 5.51 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 10.6 (264, 252) 262, 255, 250 Phenglutarimide- 258nm 5.78 264, 252 Cyclopentolate-HCI 258nm 9.91 (264, 253) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 5.64 (264, 252) 230 Fenproporex-HCI 258nm 4.36 (264, 252) Ipratropium Brom 258nm 4.36 (264, 252) Oxybutynin-HCI 258nm 5.38	HCI leate drine-HCI e-HCI HCI Bromide	1902 1227 324 3048 1530 3051 3052 3053 1910 032 1455 30104
258nm 695 240 ∠ ∠ Chlorhexidine-2H 258nm 370 236 ↑ ↑ Dimetindene male 258nm 7.25 (264, 253) 263, 255, 228 Mefenorex-HCI 258nm 18.4 (264, 252, 248) (257, 244) Mesuximide 258nm 9.82 (264, 252) 262, 254, 229 Norpseudoephed 258nm 5.51 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 10.6 (264, 252) 262, 255, 250 Phenglutarimide- 258nm 5.78 264, 252 262, 255, 250 Phenglutarimide- 258nm 9.91 (264, 252) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 4.36 (264, 252) Ipratropium Brom 258nm 4.36 (264, 252) Oxybutynin-HCI 258nm 5.38 (264, 252) Oxybutynin-HCI	HCI leate drine-HCI e-HCI HCI Bromide	1227 324 3048 1530 3051 3052 3053 1910 032 1455 30104
258nm 370 236 ↑ ↑ Dimetindene male 258nm 7.25 (264, 253) 263, 255, 228 Mefenorex-HCI 258nm 18.4 (264, 252, 248) (257, 244) Mesuximide 258nm 9.82 (264, 252) 262, 254, 229 Norpseudoephed 258nm 5.51 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 10.6 (264, 252) 262, 255, 250 Phenglutarimide- 258nm 5.78 264, 252 Cyclopentolate-H 258nm 9.91 (264, 253) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 1400 229 Mesulfen 258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.38 (264, 252) R	drine-HCI e-HCI HCI Bromide	324 3048 1530 3051 3052 3053 1910 032 1455 30104
258nm 7.25 (264, 253) 263, 255, 228 Mefenorex-HCI 258nm 18.4 (264, 252, 248) (257, 244) Mesuximide 258nm 9.82 (264, 252) 262, 254, 229 Norpseudoephed 258nm 5.51 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 10.6 (264, 252) 262, 255, 250 Phenglutarimide- 258nm 5.78 264, 252 Cyclopentolate-H 258nm 9.91 (264, 253) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.0 (267, 252) Ramipril	drine-HCI B-HCI HCI SI Bromide	3048 1530 3051 3052 3053 1910 032 1455 30104
258nm 18.4 (264, 252, 248) (257, 244) Mesuximide 258nm 9.82 (264, 252) 262, 254, 229 Norpseudoephed 258nm 5.51 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 10.6 (264, 252) 262, 255, 250 Phenglutarimide- 258nm 5.78 264, 252 Cyclopentolate-H 258nm 9.91 (264, 253) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.0 (267, 252) Ramipril	e-HCI HCI HCI Bromide	1530 3051 3052 3053 1910 032 1455 30104
258nm 9.82 (264, 252) 262, 254, 229 Norpseudoephed 258nm 5.51 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 10.6 (264, 252) 262, 255, 250 Phenglutarimide- 258nm 5.78 264, 252 Cyclopentolate-H 258nm 9.91 (264, 253) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.0 (267, 252) Ramipril	e-HCI HCI HCI Bromide	3051 3052 3053 1910 032 1455 30104
258nm 5.51 (264, 252) 262, 255, 239 Pentorex tartrate 258nm 10.6 (264, 252) 262, 255, 250 Phenglutarimide- 258nm 5.78 264, 252 Cyclopentolate-H 258nm 9.91 (264, 253) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.0 (267, 252) Ramipril	e-HCI HCI HCI Bromide	3052 3053 1910 032 1455 30104
258nm 10.6 (264, 252) 262, 255, 250 Phenglutarimide-Cyclopentolate-H 258nm 5.78 264, 252 Cyclopentolate-H 258nm 9.91 (264, 253) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 1400 229 Mesulfen 258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.0 (267, 252) Ramipril	-HCI HCI SI 3 Bromide	3053 1910 032 1455 30104
258nm 5.78 264, 252 Cyclopentolate-H 258nm 9.91 (264, 253) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 1400 229 Mesulfen 258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.0 (267, 252) Ramipril	HCI 3	1910 032 1455 30104
258nm 9.91 (264, 253) 249 Doxapram-HCI 258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 1400 229 Mesulfen 258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.0 (267, 252) Ramipril	Bromide	032 1455 30104
258nm 16.8 238 Fendiline-HCI 258nm 8.0 (263, 252) 230 Fenproporex-HCI 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 1400 229 Mesulfen 258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.0 (267, 252) Ramipril	Bromide	1455 30104
258nm 8.0 (263, 252) 230 Fenproporex-HCl 258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 1400 229 Mesulfen 258nm 5.38 (264, 252) Oxybutynin-HCl 258nm 5.0 (267, 252) Ramipril	Bromide	30104
258nm 5.64 (264, 252) Glycopyrronium E 258nm 4.36 (264, 252) Ipratropium Brom 258nm 1400 229 Mesulfen 258nm 5.38 (264, 252) Oxybutynin-HCl 258nm 5.0 (267, 252) Ramipril	Bromide	
258nm 4.36 (264, 252) Ipratropium Brom 258nm 1400 229 Mesulfen 258nm 5.38 (264, 252) Oxybutynin-HCl 258nm 5.0 (267, 252) Ramipril		2455
258nm 1400 229 Mesulfen 258nm 5.38 (264, 252) Oxybutynin-HCl 258nm 5.0 (267, 252) Ramipril	nide	2700
258nm 5.38 (264, 252) Oxybutynin-HCI 258nm 5.0 (267, 252) Ramipril		2457
258nm 5.0 (267, 252) Ramipril	3	30106
		2466
258nm 5 92 (263, 252) 230 Tribevulphenidul-		628
	-HCI 3	30120
258nm 12.8 (268, 264, 251) 261, 254, 249 Pipoxolan-HCl		2477
258nm 9.6 (262, 252) 262, 255, 249 Fenpiverinium Br	romide	2404
257nm 5.6 (263, 251) 262, 255, 248 Atropine sulfate		1901
257nm 5.0 263, 252 262, 254, 248 Atropine Methobr	romide	1906
257nm 5.2 263, 252 262, 254, 248 Atropine Methoni	itrate	1907
257nm 538 236		2482
257nm 4.3 263, 252 262, 254, 248 Butylscopolamini Bromide	ium	2435
257nm 220 247 ↓ Fominoben-HCI		023
257nm 446 242 Melitracen-HCI		2175
257nm 367 242 Melitracen mesila	ate	2176
257nm 4.6 263, 252 261, 254, 248 Oxitropium Brom	iide	2447
257nm 518 236 Tolciclate		9102
257nm 6.0 (263, 251) 230 Procyclidine-HCI		2412
257nm 3.2 (263, 251) 262, 255, 248 (↑) Scopolamine-HB	Br	1903
257nm 7.2 (263, 251) 261, 254, 248 Pethidine-HCl		113
257nm 10.5 (263, 251) 223 Propoxyphene-H	ICI	116
257nm 6.3 (263, 251) 248 Tilidine-HCl		120
257nm 420 <u>246</u> 425 259, 213 Propyphenazone	;	117
257nm 6.2 (266, 263), 251 238 Phendimetrazine		2149
257nm 7.9 Azidocillin sodiun		473
257nm 4.67 (263, 251) Methylscopolamii Bromide	inium	2459
257nm 4.91 (263, 251) Methylscopolamii Nitrate	inium	2460
257nm 9.4 (263, 251) 229 Phenylpropanolal	ımine-HCl	2731
256nm 13 243 Pentaerythrityl Te		1442
256nm 827 224 7 Propylparaben so		1208a
256nm 80 241 Tropicamide		1909
256nm 760 231 ↓ ↓ Opipramol-2HCI		2119
256nm 18,2 (262, 250) (254, 247) ↓ Glutethimide		1504
256nm 928 224 7 Propylparaben		1208
256nm 1063 224 7 Methylparaben		1207
256nm 1010 224 7 Ethylparaben		1206

Max. in methanol 256nm 10.8 (262, 225) 225 225 7 Buylparaben 132 256nm 880 220 7 Prednimustne 9122 256nm 485 220 7 Prednimustne 9122 256nm 86 226 7 Prednimustne 9122 256nm 8.6 261, 249 230 7 Prednimustne 9122 256nm 8.6 226 3 3 3 226 3 226 3 3 3 3 226 3 226 3 3	Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1% 1cm	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
ZeSemm 10.8 (262.225) 225 7						acidic	alkal.		
2256 mm		10.8	(262, 225)		225			Ephedrine-HCI	003
Z55mm			. ,		225		7	•	
Z55mm 8.6 (Z61, 249) 230 Phenmetrazin-HCI 007 285nm 872 220 J K Arsanitizacid 905 255nm 672 231 → Ketoprofen 716 255nm 545 228 N Lofepramine-HCI 214 255nm 404 221 J V Valacidovir 981 254nm 611 221 J J Acicolovir 981 254nm 300 228 Biforazole 984 254nm 300 228 Strychnine sulfate 014 254nm 300 233 J Feritionazole 291 252nm 300 239 J Feritionazole intrate 913 252nm 302 (276) 219 X 257 Methyamine phenazone 2910 252nm 302 (276) 219 X 257 Methyamine phenazone 2910 252nm 302 (276)	256nm	485			220			• •	
Z55mm 8.6 (Z61, 249) 230 Phenmetrazin-HCI 007 285nm 872 220 J K Arsanitizacid 905 255nm 672 231 → Ketoprofen 716 255nm 545 228 N Lofepramine-HCI 214 255nm 404 221 J V Valacidovir 981 254nm 611 221 J J Acicolovir 981 254nm 300 228 Biforazole 984 254nm 300 228 Strychnine sulfate 014 254nm 300 233 J Feritionazole 291 252nm 300 239 J Feritionazole intrate 913 252nm 302 (276) 219 X 257 Methyamine phenazone 2910 252nm 302 (276) 219 X 257 Methyamine phenazone 2910 252nm 302 (276)	255nm	591			240		И	Diflunisal	
255nm 872 220 ↓ k Arsamilic acid 905 255nm 672 231 → Ketprofen 716 255nm 545 228 Y Lofepramine-HCI 2145 256nm 404 221 ↓ Å Valacidovir-HCI 918 254nm 778 228 Bifonazole 984 254nm 300 228 Strychine suffate 014 254nm 374 218 Clostebol 17}-actate 228 254nm 374 218 Clostebol 17}-actate 228 256nm 300 239 ↓ Ferticonazole firate 9137 252nm 398 224 Keffilm Felbinac 747 252nm 392 (276) 219 X 257 Methylamino phenazone 290 252nm 392 (276) 219 X 257 Methylamino phenazone 291 252nm 392 (276) 219 X 257 Methylamino phenazone 291			(261, 249)						
255nm	255nm	872			220	\downarrow	Ľ	Arsanilic acid	905
255mm 545	255nm	672			231	\rightarrow	\rightarrow	Ketoprofen	716
255mm 404 221	255nm	545			228			•	
254nm 611 221 ↓ ↓ Aciclovir 981 254nm 778 228 Bl6nazole 984 254nm 300 228 Strychnine sulfate 014 254nm 374 218 Clostebol 17β-acetate 2254 254nm 370 234 Natfifne-HCI 9118 253nm 300 239 ↓ Ferticonazole nitrate 9137 252nm 998 224 Felbinac 747 252nm 998 224 Felbinac 747 252nm 392 (276) 219 X257 Methylamino phenazone 2910 252nm 1088 221 > X57 Methylamino phenazone 2910 252nm 2335 215 > Sorbic acid 1205 252nm 172 215 > Sorbic acid 1201 252nm 230 Ceflazidime 477 252nm 290 230 L Fluvoxamine maleate 2183 251nm	255nm				221	\downarrow	V	•	
254nm 778 228 Bifonazole 984 254nm 300 228 Strychnine sulfate 014 254nm 374 218 Clostebol 17β-acetate 2254 254nm 675 (291, 281) 234 Naffifine-HCl 9118 253nm 300 239 ↓ Fenticonazole nitrate 9137 252nm 998 224 Kericonazole nitrate 9137 252nm 1088 221 7 4-Hydroxybenzole acid 1205 252nm 1088 221 7 4-Hydroxybenzole acid 1205 252nm 1088 221 7 4-Hydroxybenzole acid 1205 252nm 1089 2215 7 Carbenokolone disodium salt 252nm 172 215 7 Carbenokolone disodium salt 252nm 395 230 L Kellowoximie maleate 283 252nm 396 230 L Kellowoximie maleate 283	254nm	611			221	\downarrow	\downarrow	Aciclovir	981
254nm 300 228 Strychnine sulfate 014 254nm 374 218 Clostebol 17β-acetate 2294 254nm 675 (291, 281) 234 Naftline-HCI 918 255nm 300 239 ↓ Fentionaczole nitrate 9137 252nm 392 (276) 219 X 257 Methylamino phenazone 2910 252nm 1068 221 7 4-Hydroxybenzoic acid 1211 252nm 1068 221 7 4-Hydroxybenzoic acid 1211 252nm 1083 215 36 215 36 Carbenzoione disodium 1607 252nm 395 230 Coffazidime 477 252nm 290 230 K Fluvoxamine maleate 218 251nm 396 230 230 K Fluvoxamine maleate 218 221 221 47 220 221 222 223 222 223 222 223 222	254nm	778			228			Bifonazole	
254nm 374 218 Clostebol 17β-acetate 2254 254nm 675 (291, 281) 234 Natifirie-HCI 9118 253nm 300 239 ↓ Fenticonazole nitrate 9137 252nm 998 224 Febbinac 747 252nm 1088 221 ¬ 4-Hydroxybenzoic acid 1205 252nm 2335 215 → Sorbic acid 1211 252nm 2335 215 → Sorbic acid 1211 252nm 2335 230 Caffacidime 477 252nm 395 230 Caffacidime 477 252nm 290 230 K K Fluvoxamine maleate 2183 251nm 398 (288, 279) 234 Carteoloi-HCI 2027 251nm 308 (288, 279) 234 Carteoloi-HCI 2027 251nm 309 (280) 220 X K 4-Aminophenazone 2912 29117 29117 29					+			Strychnine sulfate	
255nm 675 (291, 281) 234 Naftifine-HCI 9118 253nm 300 239 ↓ Fenticonazole nitrate 9137 252nm 998 224 Felbinac 747 252nm 392 (276) 219 X 257 Methylamino phenazone 2910 252nm 1068 221 ¬ 4 - Hydroxybenzole acid 1201 252nm 2335 215 → Sorbic acid 1211 252nm 172 215 → Gorbenoxolone disodium salt 1607 252nm 395 230 ∠ Le Fluvoxamine maleate 2183 251nm 398 (288, 279) 234 Carteolol-HCI 2027 251nm 398 (288, 279) 234 Carteolol-HCI 2027 251nm 230 Le Fluvoxamine maleate 2183 251nm 238 237 → Clomipramine-HCI 2104 251nm 230 Le Fluvoxamine maleate 2104 251nm 1300 220 X					218			•	
253nm 300 239	254nm	675	(291, 281)		234			•	
252nm 998 224 X 57 Methylamino phenazone 2910 252nm 392 (276) 219 X 257 Methylamino phenazone 2910 252nm 1068 221 A 4-Hydroxybenzoic acid 1205 252nm 2335 215 → Sorbic acid 1211 252nm 395 230 Cafbacidine 4.77 252nm 395 230 Cefazidine 4.77 252nm 395 230 Cafteoloi-HCl 2027 251nm 308 (288, 279) 234 Carteoloi-HCl 2027 251nm 274 230 Imipramine-HCl 2104 251nm 234 237 → Clomipramine-HCl 2104 251nm 549 (280) 220 X X 4-Arminophenazone 2912 251nm 549 (280) 220 X X 4-Arminophenazone 2912 251nm 1300 223 235 222 Ciondifiazione<	253nm		(- , - ,			1			
252nm 392 (276) 219 X 257 Methylamino phenazone 2910 252nm 1068 221 7 4-Hydroxyberacic acid 1205 252nm 2335 215 → Sorbic acid 1211 252nm 172 215 → Cardenoxolone disodium salt 1807 252nm 395 230 Ceffazidime 477 252nm 290 230 L Fluvoxamine maleate 2183 251nm 308 (288, 279) 234 Carteolol-HCI 2027 251nm 238 237 → Clomipramine-HCI 2117 251nm 238 237 → Clomipramine-HCI 2117 251nm 1300 220 X L 4-Aminophenazone 2912 251nm 562 293, 284, 228 235, 222 Cinnarizine 016 250nm 176 224 7 Clomethiazole edisilate 216 250nm 176 224<									
252nm			(276)			X 257			
252nm			(=)			1	7	•	
252nm 172	-					\rightarrow		•	
252nm 290 230 ∠ ∠ Fluvoxamine maleate 2183 251nm 308 (288, 279) 234 Carteolol-HCI 2027 251nm 274 230 Imipramine-HCI 2104 251nm 238 237 → Clomipramine-HCI 2117 251nm 549 (280) 220 X ∠ 4-Aminophenazone 2912 251nm 1300 223 ↓ ¾ Azapropazone dihydrate 715 251nm 562 293, 284, 228 235, 222 Cinnarizine 0 16 250nm 175 224 ౫ ↓ Clomethiazole edisilate 2160 250nm 296 230 Desipramine-HCI 2118 250nm 410 232 Noxiptiine-HCI 2147 249nm 740 217 ∠ ∠ p-Latophenetide 142 249nm 566 230 ¸ Allopurinol 2501 249nm 566 230<							\rightarrow	Carbenoxolone disodium	
251nm 308 (288, 279) 234 Carteolol-HCl 2027 251nm 274 230 Impramine-HCl 2104 251nm 238 237 → Clomipramine-HCl 2117 251nm 549 (280) 220 X	252nm	395			230			Ceftazidime	477
251nm 308 (288, 279) 234 Carteolol-HCl 2027 251nm 274 230 Impramine-HCl 2104 251nm 238 237 → Clomipramine-HCl 2117 251nm 549 (280) 220 X	252nm	290			230	L	Ľ	Fluvoxamine maleate	
251nm 274 230 Imipramine-HCI 2104 251nm 238 237 → Clomipramine-HCI 2117 251nm 549 (280) 220 X X 4-Aminophenazone 2912 251nm 1300 223 J J Azapropazone dihydrate 715 251nm 562 293, 284, 228 235, 222 Cinnarizine 016 250nm 175 224 A J Clomethiazole edisilate 2160 250nm 296 230 Desipramine-HCI 2118 250nm 410 232 Noxiptiline-HCI 2148 249nm 740 217 ∠ ∠ p-Lactophenetide 142 249nm 563 219 → ∠ Amfepramone-HCI 3062 249nm 566 230 J Allopurinol 2501 249nm 566 230 J Allopurinol 2501 249nm 864 216 <			(288, 279)						
251nm 238 237 → Clomipramine-HCl 2117 251nm 549 (280) 220 X ∠ 4-Aminophenazone 2912 251nm 1300 223 ↓ Napropazone dihydrate 715 251nm 562 293, 284, 228 235, 222 Cinnarizine 016 250nm 175 224 ७ Ciomethiazole edisilate 2160 250nm 296 230 Desipramine-HCl 2118 250nm 410 232 Noxiptiline-HCl 2118 250nm 740 217 ∠ ∠ P-Lactophenetide 142 249nm 740 217 ∠ ∠ P-Lactophenetide 142 249nm 238 235 ↓ Trimipramine maleate 2136 249nm 238 235 ↓ Trimipramine maleate 2136 249nm 566 230 Nallopurinol 2501 249nm 864 216 ↓ Phenacetin 114 249nm 8.3 (222) 238, 217 K	+		. ,						
251nm 1300 223 ↓ N Azapropazone dihydrate 715 251nm 562 293, 284, 228 235, 222 Cinnarizine 016 250nm 175 224 7 ↓ Clomethiazole edisilate 2160 250nm 296 230 Desipramine-HCI 2118 250nm 410 232 Noxiptiline-HCl 2147 249nm 740 217 ∠ ∠ p-Lactophenetide 142 249nm 563 219 → ∠ Amfepramone-HCl 3062 249nm 238 235 ↓ Trimipramine maleate 2136 249nm 566 230 µ Allopurinol 2501 249nm 709 216 ↓ Bucetin 104 249nm 8.3 (222) 238,217 下 Ethosuximide 1528 249nm 8.3 (222) 238,217 下 Ethosuximide 1528 249nm <t< td=""><td>251nm</td><td>238</td><td></td><td></td><td>237</td><td></td><td>\rightarrow</td><td>•</td><td>2117</td></t<>	251nm	238			237		\rightarrow	•	2117
251nm 562 293, 284, 228 235, 222 Cinnarizine 016 250nm 175 224 7 ↓ Clomethiazole edisilate 2160 250nm 296 230 Desipramine-HCI 2118 250nm 410 232 Noxiptiline-HCI 2147 249nm 740 217 ∠ ∠ p-Lactophenetide 142 249nm 563 219 → ∠ Amfepramone-HCI 3062 249nm 238 235 ↓ Trimipramine maleate 2136 249nm 566 230 ⅓ Allopurinol 2501 249nm 864 216 ↓ ↓ Bucetin 104 249nm 8.3 (222) 238, 217 K Ethosuximide 1528 249nm 624 220 Metamfepramone-HCI 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ Probenecid	251nm	549	(280)		220	Х	Ľ	4-Aminophenazone	2912
250nm 175 224 7 ↓ Clomethiazole edisilate 2160 250nm 296 230 Desipramine-HCI 2118 250nm 410 232 Noxiptiline-HCI 2147 249nm 740 217 ⊭ ⊭ p-Lactophenetide 142 249nm 563 219 → ⊭ Amfepramone-HCI 3062 249nm 563 230 ⇒ Allopurinol 2501 249nm 566 230 ⇒ Allopurinol 2501 249nm 709 216 ⇒ ⇒ Bucetin 104 249nm 864 216 ⇒ ⇒ Phenacetin 114 249nm 8.3 (222) 238, 217 ⋈ Ethosuximide 1528 249nm 624 220 Metamfepramone-HCI 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ⇒ Practolol 2006<	251nm	1300			223	\downarrow	Z	Azapropazone dihydrate	715
250nm 175 224 7 ↓ Clomethiazole edisilate 2160 250nm 296 230 Desipramine-HCI 2118 250nm 410 232 Noxiptiline-HCI 2147 249nm 740 217 ⊭ ⊭ p-Lactophenetide 142 249nm 563 219 → ⊭ Amfepramone-HCI 3062 249nm 566 230 ⊭ Allopurinol 2501 249nm 566 230 ⊭ Allopurinol 2501 249nm 709 216 ⊭ Bucetin 104 249nm 864 216 ⊭ Phenacetin 114 249nm 8.3 (222) 238, 217 ⋈ Ethosuximide 1528 249nm 624 220 Metamfepramone-HCI 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ⊭ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530	251nm	562	293, 284, 228		235, 222			Cinnarizine	016
250nm 296 230 Desipramine-HCI 2118 250nm 410 232 Noxiptiline-HCI 2147 249nm 740 217 ∠ ∠ p-Lactophenetide 142 249nm 563 219 → ∠ Amfepramone-HCI 3062 249nm 566 235 ↓ Trimipramine maleate 2136 249nm 566 230 ¾ Allopurinol 2501 249nm 709 216 ↓ ↓ Bucetin 104 249nm 864 216 ↓ ↓ Phenacetin 114 249nm 8.3 (222) 238,217 ⋉ Ethosuximide 1528 249nm 624 220 Metamfepramone-HCI 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) <td< td=""><td>250nm</td><td>175</td><td></td><td></td><td>224</td><td>7</td><td>\downarrow</td><td>Clomethiazole edisilate</td><td>2160</td></td<>	250nm	175			224	7	\downarrow	Clomethiazole edisilate	2160
250nm 410 232 Noxiptiline-HCI 2147 249nm 740 217 ∠ ∠ p-Lactophenetide 142 249nm 563 219 → ∠ Amfepramone-HCI 3062 249nm 238 235 ↓ Trimipramine maleate 2136 249nm 566 230 ⅓ Allopurinol 2501 249nm 709 216 ↓ ↓ Bucetin 104 249nm 864 216 ↓ ↓ Phenacetin 114 249nm 8.3 (222) 238, 217 ⋉ Ethosuximide 1528 249nm 624 220 Metamfepramone-HCI 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 247nm 429 225 ∠ ∠ Atovastatin calcium salt, trihydr 1816 247nm 808	250nm							Desipramine-HCI	
249nm 740 217 ∠ ∠ p-Lactophenetide 142 249nm 563 219 → ∠ Amfepramone-HCl 3062 249nm 238 235 ↓ Trimipramine maleate 2136 249nm 566 230 IJ Allopurinol 2501 249nm 709 216 ↓ ↓ Bucetin 104 249nm 864 216 ↓ ↓ Phenacetin 114 249nm 8.3 (222) 238, 217 ⋉ Ethosuximide 1528 249nm 624 220 Metamfepramone-HCl 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ∠								•	
249nm 563 219 → ८ Amfepramone-HCI 3062 249nm 238 235 ↓ Trimipramine maleate 2136 249nm 566 230 ы Allopurinol 2501 249nm 709 216 ↓ ↓ Bucetin 104 249nm 864 216 ↓ ↓ Phenacetin 111 249nm 8.3 (222) 238, 217 下 Ethosuximide 1528 249nm 624 220 Metamfepramone-HCl 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ८ ८ Atorvastatin calcium salt, trihydr 1816 247nm 860 216 ↓ □ Paracetamol 112 247nm 2 232						Ľ	Ľ	•	
249nm 238 235 ↓ Trimipramine maleate 2136 249nm 566 230 ы Allopurinol 2501 249nm 709 216 ↓ Bucetin 104 249nm 864 216 ↓ Phenacetin 114 249nm 8.3 (222) 238, 217 K Ethosuximide 1528 249nm 624 220 Metamfepramone-HCl 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ८ ८ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 2 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm <t< td=""><td>-</td><td></td><td></td><td></td><td></td><td>\rightarrow</td><td>Ľ</td><td></td><td></td></t<>	-					\rightarrow	Ľ		
249nm 566 230 ☑ Allopurinol 2501 249nm 709 216 ☑ Bucetin 104 249nm 864 216 ☑ Phenacetin 114 249nm 8.3 (222) 238, 217 ☒ Ethosuximide 1528 249nm 624 220 Metamfepramone-HCl 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ☑ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ☒ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ☑ Flurbiprofen 724 247nm 850 216 ☑ Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 436 226 ↑ ☒ Azapetine phosphate 2032 247nm	249nm	238			235			Trimipramine maleate	
249nm 709 216 ↓ ↓ Bucetin 104 249nm 864 216 ↓ ↓ Phenacetin 114 249nm 8.3 (222) 238, 217 ⋉ Ethosuximide 1528 249nm 624 220 Metamfepramone-HCI 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ८ ८ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ ¥ Paracetamol 112 247nm 457 224 Piperylone 144 247nm 457 224 Piperylone maleate 145 247nm 436 226 ↑ ८ </td <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td>И</td> <td>•</td> <td></td>							И	•	
249nm 864 216 ↓ ↓ Phenacetin 114 249nm 8.3 (222) 238, 217 ⋉ Ethosuximide 1528 249nm 624 220 Metamfepramone-HCI 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ८ ८ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ у Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 436 226 ↑ ८ Azapetine phosphate 2032 247nm 462 224 Cyclofenil						\downarrow	\downarrow	•	
249nm 8.3 (222) 238, 217 N Ethosuximide 1528 249nm 624 220 Metamfepramone-HCl 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ∠ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ ☑ Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242	249nm	864			216		<u> </u>	Phenacetin	
249nm 624 220 Metamfepramone-HCI 30107 249nm 309 243 Valsartan 1335 248nm 610 215 ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ∠ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ J Paracetamol 112 247nm 457 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 436 222 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242			(222)		ļ	+			
249nm 309 243 Valsartan 1335 248nm 610 215 ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ८ ८ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ у Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 436 222 Piperylone maleate 145 247nm 436 226 ↑ ८ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242			()						
248nm 610 215 ↓ Practolol 2006 248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ∠ ∠ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ □ Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 290 222 Piperylone maleate 145 247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242					ļ			· · · · · · · · · · · · · · · · · · ·	
248nm 22.5 (264, 258, 252) (257, 244) Mesuximide 1530 248nm 332 227 234 Probenecid 3015 247nm 429 225 ∠ ∠ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ ✓ Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 290 222 Piperylone maleate 145 247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242						\downarrow	\downarrow		
248nm 332 227 234 Probenecid 3015 247nm 429 225 ∠ ∠ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ ✓ Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 290 222 Piperylone maleate 145 247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242			(264, 258, 252)			1	•		
247nm 429 225 ∠ Atorvastatin calcium salt, trihydr 1816 247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ ☑ Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 290 222 Piperylone maleate 145 247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242			·						
247nm 808 222 ↓ Flurbiprofen 724 247nm 850 216 ↓ ☑ Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 290 222 Piperylone maleate 145 247nm 436 226 ↑ ८ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242						Ľ	Ľ		
247nm 850 216 ↓ ☑ Paracetamol 112 247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 290 222 Piperylone maleate 145 247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242								•	
247nm Z 232 Amfetaminil 018 247nm 457 224 Piperylone 144 247nm 290 222 Piperylone maleate 145 247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242	+						И	· · · · · · · · · · · · · · · · · · ·	
247nm 457 224 Piperylone 144 247nm 290 222 Piperylone maleate 145 247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242	+					1			
247nm 290 222 Piperylone maleate 145 247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242									
247nm 436 226 ↑ ∠ Azapetine phosphate 2032 247nm 462 224 Cyclofenil 2242								· · ·	
247nm 462 224 Cyclofenil 2242	-					↑	Ľ		
						· ·			
7 -7 (CATANCIAN CANCIAN SPIN 17:14	247nm	1710				\rightarrow	\rightarrow	Sorbic acid potassium salt	1234

Distinct long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Changionis		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
246nm	104			230			Benzyl isothiocyanate	970
246nm	424			219	1		Acetohexamide	501
246nm	423	238, 230	631	244, 233,<210	<u> </u>	V	Lovastatin	1818
246nm	Z		031				Simvastatin	1815
		238, 229		243, 233, <210		N		
245nm	475			214		<u> </u>	Nimustine-HCI	9106
245nm	417			230	X	Z	Oxazolam	2627
245nm	393			221	1	Ľ	Sultiame	1542
245nm	520			<210	\rightarrow	\rightarrow	Metandienone	1015
245nm	84			230	7		Haloperidol	2112
245nm	259			(238)			Trimetozine	2608
245nm	502			225	Ľ	7	Sulfinpyrazone	2507
245nm	469			<210			Ascorbic acid	2801
245nm	455	(215)				X 260	Disulfiram	3042
245nm	267			224	71	Z	Eprazinone-2HCl	813
245nm	522	227	500	235, 219	K		Almitrine dimesilate	031
245nm	340			215			Halometasone	1052
244nm	335			227			Moperone-HCI	2168
244nm	287			220			Pipamperone-2HCl	2170
244nm	290			223	\rightarrow		Trifluperidol-HCl	2135
244nm	448			223	L L	7	Kebuzone	712
244nm	296			213			Methylprednisolone 21- (sodium succinate)	1022
244nm	412			214			Prednylidene	1023
244nm	327	(227)		214			Prifinium Bromide	2461
244nm 243nm	349	(221)		<210			Deflazacort	1055
243nm	508			218				827
				221			Dropropizine	
243nm	413						Melperone-HCI	2167
243nm 243nm	359 302			<220 <210			Iohexol Methylprednisolone 21-	3080 1037
							cipionate	1037
243nm	482			223	\downarrow	7	Phenylbutazone	701
243nm	360			210			Methylprednisolone 21- acetate	1018
243nm	370			211			Prednisolone 21-acetate	1017
243nm	428			228			Probucol	1805
243nm	347			<220			Budenoside	1040
243nm	363			<220			Desonide	1042
243nm	684			227			Flumazenil	30123
243nm	268	225	330	237	Х		Terconazole	9125
242nm	370			<210			Flunisolide	1056
242nm	297			224			Tiamenidine-HCI	617
242nm	500			221	1	Z	Aminoglutethimide	3091
242nm	416			(210)	*	_	Prednisolone	1011
2421111 242nm	404			(210)			Methylprednisolone	1011
242nm	500			<210			Hydroxyprogesterone	2224
				<220			Fluocortolone	
242nm	428							1044
242nm	317			<220			Prednicarbate	1054
241nm	490			000			Ketazolam	2621
241nm	460			222			Benfotiamine phosphate	2816
241nm	387			<210			Metenolone 17β-acetate	2236
241nm	536			<210			Ethisterone	2218
241nm	480			220	←	7	Oxyphenbutazone	710
241nm	1060			214	↓	\downarrow	Acetanilide	101
241nm	437			<210			Aldosterone	1020

Distinct long wave length	E 1% 1cm	Further maxima (nm) in methanol	E 1% 1cm	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
241nm	564			<210			Testosterone	2220
241nm	311	227, 221		234	\rightarrow		Bromperidol	2180
241nm	379			<220			Fluocortin Butyl	1043
241nm	358			<220			Hydrocortisone 21-acetate 17-propionate	1038
241nm	402			<220			Hydroxyprogesterone 17α- caproate	2255
241nm	368			<220			Paramethasone 21-acetate	1053
241nm	545			<210			Testolactone	2253
240nm	553			222	Z		Amiodarone-HCI	1445
240nm	416			<210			Medroxyprogesterone 17α-acetate	2234
240nm	381			<210			Fluprednisolone 21-acetate	1007
240nm	510			<210	\rightarrow	\rightarrow	Norethisterone 17β-acetate	2211
240nm	391			<210			Hydrocortisone 21-acetate	1009
240nm	573			<210	\rightarrow	\rightarrow	Norethisterone	2212
240nm	407			<210	<u> </u>	,	Testosterone 17β-cipionate	2209
240nm	554			<210	\rightarrow	\rightarrow	Norgestrel	2214
240nm	533			<210	7		Methyltestosterone	2206
240nm	540			<210			Progesterone	2208
240nm	384			210			Triamcinolone	1024
240nm	312			209			Triamcinolone 16α, 21- diacetate	1024
240nm	291			210			Triamcinolone hexacetonide	1027
240nm	503			226	X		Brotizolam	1544
240nm	546			<210			Norgestimate	2250
239nm	537			223			Diltiazem-HCl	1437
239nm	296			210			Beclometasone 17α, 21-	1437
239nm	747			215			dipropionate Benorilate	1030 731
								731
239nm	259			<210			Betamethasone 21- (disodium phosphate) Dexamethasone 21-(3,6,9-	1031
239nm	268			<210			trioxa)-undecanoate	1033
239nm	865			220		L	Moroxydine-HCl	996
239nm	402			<210			Nandrolone 17β-decanoate	2237
239nm	432			218			Nandrolone 17β- phenpropionate	2238
239nm	472			228			Nortriptyline-HCI	2127
239nm	410			(210)			Fluorometholone	1013
239nm	441			<210			Hydrocortisone	1008
239nm	444			<210			Desoxycortone 21-acetate	1004
239nm	386			208			Dexamethasone	1019
239nm	328			210			Flumetasone 21-pivalate	1021
239nm	380			209			Prednisone 21-acetate	1016
239nm	347			208			Triamcinolone 16α , 17α -acetonide	1025
239nm	490			220	\	-	Etomidate	1545
239nm	569			<210			Gestodene	2249
239nm	423			<210			Gestonorone caproate	2245
238nm	219			222			Benazepril-HCl	629
238nm	339			<210			Clobetasol 17α-propionate	1032
238nm	456			<210			Fludrocortisone	1034
238nm	411			<210			Fludrocortisone 21-acetate	1035

Max. in methanol 229	Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1%	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
238nm	Max. in					acidic	alkal.		
238nm 365		441			229	\downarrow	\downarrow	Amitriptylinoxide	2158
238nm 387 212 Betamethassone 1001 238nm 429 <210	238nm	343			<210			Diflorasone 17α, 21-	
228nm 422 < 210	238nm	365			<210			Dexamethasone 21-acetate	1005
238mm 429	238nm	387			212			Betamethasone	1001
238mm	238nm	422			<210			Fluprednisolone	1006
238mm	238nm	429			<210			Desoximetasone	1003
238nm	238nm	421			(210)			Prednisone	1012
238nm 505 (210) 7 N ω-Santonin 204 238nm 308 <220 Amcinonide 1039 238nm 499 <210 Fluoxymesterone 2248 238nm 400 <220 Fluorymesterone 2248 238nm 361 <220 Fluorymesterone 1049 238nm 367 <220 Fluorednidene 21-acetate 1050 238nm 377 <220 Halcinonide 1051 237nm 375 221 ↓ 7 Suxibuzone 729 237nm 509 219 ∠ ∠ Bumadizone calcium salt 722 237nm 315 267 180 260, 216 ↓ Z Cefoxtin sodium salt 439 237nm 328 <220 Clocortolone 21-pivalate 1041 237nm 395	238nm	471			<210		\rightarrow X	Spironolactone	2502
238nm 308	238nm	464			227		\downarrow	Amitriptyline-HCI	2101
238nm 499	238nm	505			(210)	7	K	α-Santonin	204
238mm 361	238nm	308			<220			Amcinonide	1039
238nm 361	238nm	499			<210			Fluoxymesterone	2248
238nm 377	238nm	400			<220			Fluprednidene	1049
237nm 375 221 ↓ β Suxibuzone 729 237nm 509 219 ⊭ ⊭ Bumadizone calcium salt 722 237nm 315 267 180 260, 216 ↓ Z Cefoxitis sodium salt 438 237nm 328 <220	238nm	361			<220			Fluprednidene 21-acetate	1050
237nm 509 219 ∠ ∠ Bumadizone calcium salt 722 237nm 315 267 180 260,216 ↓ Z Cefoxitin sodium salt 439 237nm 328 <<220	238nm	377			<220			Halcinonide	1051
237nm 315 267 180 260, 216 ↓ Z Cefoxitin sodium salt 439 237nm 328 <220	237nm	375			221	\downarrow	7	Suxibuzone	729
237nm 328	237nm	509			219	Ľ	L	Bumadizone calcium salt	722
237nm 395 Cortisone 21-acetate 1002 237nm 347 <210	237nm	315	267	180	260, 216	\downarrow	Z	Cefoxitin sodium salt	439
237nm 395 Cortisone 21-acetate 1002 237nm 347 <210	237nm	328			<220			Clocortolone 21-pivalate	1041
237nm 358 <220	237nm	395							1002
237nm 335	237nm	347			<210			Diflucortolone 21-valerate	1048
237nm 372 <220	237nm	358			<220			Fluocinolone acetonide	1045
237nm 510 (245, 231) <210	237nm	335			<220				1046
236nm 350	237nm	372			<220			Fluticasone 17-propionate	1057
236nm 345 227 № ↑ Fosfestrol 2233 236nm 330 <210	237nm	510	(245, 231)		<210	\downarrow	\downarrow	Pravastatin sodium salt	1814
236nm 330 <210	236nm	350			<210			Fludroxycortide	1036
236nm 908 218 ↓	236nm	345			227	И	↑	Fosfestrol	2233
236nm 409 226 ♂ Diethylstilbestrol 2201 235nm 397 215 → Cefotaxime sodium salt 424 235nm 330 215 → Cefpodoxime proxetil 498 235nm Z 219 → Z Ceftizoxime sodium salt 454 235nm 212 224 ↓ Z Dapiprazole-HCI 2033 234nm 350 ¬ Z Cefmenoxime-HCI 461 233nm 374 (212) ¾ Cephaloridine 404 233nm 220 Tiemonium Mesilate 1912 232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210	236nm	330			<210			Clobetasone 17α-butyrate	1028
235nm 397 215 → Cefotaxime sodium salt 424 235nm 330 215 → Cefpodoxime proxetil 498 235nm Z 219 → Z Ceftizoxime sodium salt 454 235nm 212 224 ↓ Z Dapiprazole-HCI 2033 234nm 350 ¬ Z Cefmenoxime-HCI 461 233nm 374 (212) ¾ Cephaloridine 404 233nm 220 Timenoium Mesilate 1912 232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210	236nm	908			218	\downarrow	Ľ	Metformin-HCI	513
235nm 397 215 → Cefotaxime sodium salt 424 235nm 330 215 → Cefpodoxime proxetil 498 235nm Z 219 → Z Ceftizoxime sodium salt 454 235nm 212 224 ↓ Z Dapiprazole-HCI 2033 234nm 350 ¬ Z Cefmenoxime-HCI 461 233nm 374 (212) ¾ Cephaloridine 404 233nm 374 (212) ¾ Cephaloridine 404 233nm 320 220 Tiracrillin disodium salt 470 232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210	236nm	409			226		7	Diethylstilbestrol	2201
235nm 330 215 → Cefpodoxime proxetil 498 235nm Z 219 → Z Ceftizoxime sodium salt 454 235nm 212 224 ↓ Z Dapiprazole-HCI 2033 234nm 350 ¬ Z Cefmenoxime-HCI 461 233nm 374 (212) ¬ Cephaloridine 404 233nm 220 220 Tiemonium Mesilate 1912 232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210	235nm	397			215	\rightarrow		·	424
235nm Z 219 → Z Ceftizoxime sodium salt 454 235nm 212 224 ↓ Z Dapiprazole-HCI 2033 234nm 350 ¬ Z Cefmenoxime-HCI 461 233nm 374 (212) ¾ Cephaloridine 404 233nm 220 Tiemonium Mesilate 1912 232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210	235nm	330			215			Cefpodoxime proxetil	498
235nm 212 224 ↓ Z Dapiprazole-HCI 2033 234nm 350 7 Z Cefmenoxime-HCI 461 233nm 374 (212) ¥ Cephaloridine 404 233nm 220 220 Tiemonium Mesilate 1912 232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210	235nm	Z			219		Z	· ·	
234nm 350 7 Z Cefmenoxime-HCl 461 233nm 374 (212) □ Cephaloridine 404 233nm 220 220 Tiemonium Mesilate 1912 232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210									
233nm 374 (212) ☑ Cephaloridine 404 233nm 220 Tiemonium Mesilate 1912 232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210	-					+			
233nm 220 Tiemonium Mesilate 1912 232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210					(212)				404
232nm 143 (228) Z Ticarcillin disodium salt 470 232nm 602 (215) → Clioxanide 208 231nm 278 <210					` '				1912
232nm 602 (215) → Clioxanide 208 231nm 278 <210						Z			470
231nm 278 <210	ļ ————————————————————————————————————				+ · · · · · · · · · · · · · · · · · · ·		\rightarrow		208
231nm Z 215 ↓ Y Procarbazine-HCI 960 231nm 575 (276, 265) 212 ↓ Chlorpropamide 504 231nm 412 226 Pizotifen maleate 341 231nm 330 <210									9132
231nm 575 (276, 265) 212 ↓ Chlorpropamide 504 231nm 412 226 Pizotifen maleate 341 231nm 330 <210						\downarrow			960
231nm 412 226 Pizotifen maleate 341 231nm 330 <210			(276, 265)						504
231nm 330 <210			,,				•		
230nm 603 (266) 216 Flupentixol-2HCl 2163 230nm 578 215 ∠ Bezafibrate 1803 229nm 498 (274) 211 Glisoxepide 515 229nm 260 219 ∠ X Prilocaine-HCl 1719 229nm 442 (259) 217 ↓ Chlorcyclizine-HCl 347 228nm 543 <									468
230nm 578 215 ∠ Bezafibrate 1803 229nm 498 (274) 211 Glisoxepide 515 229nm 260 219 ∠ X Prilocaine-HCl 1719 229nm 442 (259) 217 ↓ ↓ Chlorcyclizine-HCl 347 228nm 543 <			(266)					•	2163
229nm 498 (274) 211 Glisoxepide 515 229nm 260 219 ∠ X Prilocaine-HCI 1719 229nm 442 (259) 217 ↓ Chlorcyclizine-HCI 347 228nm 543 <200			, ,			ĸ		•	1803
229nm 260 219 ∠ X Prilocaine-HCl 1719 229nm 442 (259) 217 ↓ Chlorcyclizine-HCl 347 228nm 543 <200			(274)						515
229nm 442 (259) 217 ↓ ↓ Chlorcyclizine-HCl 347 228nm 543 <200			,			Ľ	Х		1719
228nm 543 <200			(259)						347
228nm Z Z Z Cefoperazone sodium salt 451			, ,			+	•	·	3114
						+	Z		
228nm 683 272 43 215 Benzyl benzoate 30147	228nm	683	272	43	215	1	=	'	30147

Distinct long wave length	E 1%	Further maxima (nm) in methanol	E 1% 1cm	Minima (nm) in methanol	Chang		Pharmaceutical substance	Spectr. No.
Max. in methanol					acidic	alkal.		
228nm	403			214	\downarrow	Х	Paclitaxel	9140
228nm	591			212			Butalamine-HCl	3096
227nm	629			213			Glimepiride	517
227nm	4.7			213			Lactitol monohydrate	1623
227nm	9.0			216			Trimethadione	1543
226nm	2400			<210	7		Altretamine	982
225nm	840			<220			Dibenzepin-HCI	2122
224nm	960			210			Halazepam	2625
224nm	146			<210	7		Stanozolol	2252
222nm	624	(268)			7		Benzoic acid sodium salt	2902a
221nm	1200			211	Z		Triazolam	1539
221nm	179						Metildigoxin	1408
221nm	1290			<210			Alprazolam	2624
221nm	1300	(273)		<220	И		Almitrine	030
221nm	86						Auranofin	739
221nm	492			<220			Tiemonium Iodide	1911
220nm	352				7		Tosylchloramide sodium	1213
220nm	431	(259)		216			Buclizine-2HCl	346
219nm	186			<200			Acetyldigoxin, α-isomer	1425
219nm	186			<200			Acetyldigoxin, β-isomer	1426
219nm	810			<210	\downarrow		Cimetidine	1608
219nm	150			<200			Lanatoside C	1429
219nm	306						Mupirocin	484
218nm	238					7	Methohexital	1540
218nm							Midazolam	1541
218nm	226						Digoxin	1405
218nm	280			<200			Peruvoside	1432
218nm	153						Deslanoside	1459
217nm	206						Strophanthin-G	1467
216nm	5.4			<200			Dicycloverine-HCI	2438
216nm	195						Digitoxin	1404
216nm	243						Pilocarpine-HCl	3029

Pharmaceutical substances showing no distinct <u>maximum</u> of absorption

A	540	litin .	20457
Acarbose Acecarbromal	516 2616	Levocarnitine	30157 425
	3092	Lincomycin hydrochloride	635
Acetylcysteine Aldrin	3027	Losartan potassium salt Lynestrenol	2221
Allylestrenol	2228	Aspartic acid magnesium salt	3083
Amikacin	472	Magnesium gluconate	3115
Arginine	30133	Magnesium glutamate hydrobromide	3082
Arginine hydrochloride	30134	Malathion	3084
Asparagin	30141	Memantine hydrochloride	2443
Aspartic acid	30140	Mesna	30159
Aspartic acid iron salt	30155	Methenamine	1224
Aurothioglucose	740	Methionin, L-	30135
Azithromycin dihydrate	4103	Methyldemeton	3033
Bencyclane fumarate	020	Methylpentynol	2622
Biotin	2817	Metrifonate	225
Bromisoval	1525	Netilmicin sulfate	466
Bromociclen	3031	Oxiconazole nitrate	9121
Calcium Gluconate	334	Pamidronic acid disodium salt	2514
Calcium Pantothenate	2815	Pancuronium Bromide	2473
Camphen	2472	Penicillamine	3088
Captopril	612	Perhexiline maleate	1457
Carbocisteine	818	Perindopril erbumine	632
Carbromal	1526	Pinen, α-isomer	2474
Cerivastatin sodium salt	1817	Pinen, β-isomer	2475
Chenodeoxycholic acid	3097	Pipecuronium Bromide	2476
Chlorobutanol	1228	Piracetam	013
Cholic acid	30168	Piperacillin sodium salt	456
Ciclosporin	3098	Piperazine	219
Eucalyptol	30153	Piperazin adipate	220
Citric acid monohydrate	3074	Polyethylene Glykol 400	3110
Clemastine fumarate	318	Polyethylene Glykol 6000	3112
Clindamycin hydrochloride	422	Propylene Glykol	3109
Colistin sulfate	480	Protirelin	3025
Crotamiton	1112	Rocuronium Bromide	2478
Deferoxamine	3024	Roxithromycin	4101
Desogestrel	2231	Sparteine sulfate pentahydrate	1443
Dexpanthenol	2820	Stearic acid	3113
Dieldrin	3026	Streptomycin sulfate	426
Dihydrostreptomycin sulfate	427	Sulbactam sodium salt	30116
Dimercaprol	3041	Suxamethonium Chloride	2425
Docusate sodium	1606	Tartaric acid	3090
Ethinamate	1524	Taurolidine	9144
Fexofenadine hydrochloride	354	Tazobactam	4102
Fosfomycin disodium salt	499	Temocillin disodium salt	9105
Fusidic acid sodium salt	485	Thiotepa	9101
Hexetidine	1223	Tranexamic acid	3089
Histamine dihydrochloride	3079	Trofosfamide	9103
Hydroxycarbamide	993	Ursodeoxycholic acid	30117
Ifosfamide	953	Valdetamide	2614
lloprost tromethamine	30169	Valproic acid sodium salt	1532
Isosorbide dinitrate	1424	Vecuronium Bromide	2481
Isosorbide mononitrate	1439	Vigabatrin	1548
Aspartic acid potassium salt	3081	Viminol	136

Pharmaceutical substances without useful UV absorption >220 nm

Acetylcysteine	3092	Isoleucine	30131
Alanine	30132	Isometheptene hydrochloride	2718
Amantadine hydrochloride	3036	Kanamycin sulfate	444
Aminocaproic acid	3093	Lactose	3107
Arginine	30133	Lactulose	30113
Arginine hydrochloride	30134	Leucine	30137
Asparagin	30141	Levoglutamide	30139
Aspartic acid	30140	Lindane	30105
Azelaic acid	1115	Lyapolate Sodium	3032
Busulfan	950	Lysine acetate	30128
Carbachol	1953	Lysine hydrochloride	30127
Carisoprodol	2453	Mannitol	3106
Cholesterol	2904	Mecamylamine hydrochloride	2442
Cholic acid	30168	Meprobamate	2615
Cyclamat-Calcium	3111	Norethynodrel	2219
Cyclamate Sodium	3117	Norleucine	30130
Cyclophosphamide monohydrate	952	Paromomycin	489
Cysteine	30142	Pentetrazol	005
Cystine	30143	Proline	30124
Dextrose	3030	Propylhexedrine hydrochloride	3054
Ethambutol dihydrochloride	949	Saccharose	3108
Framycetin sulfate	471	Serine	30136
Fruktose	3028	Sitosterol, β-	1819
Gentamicin sulfate	464	Sorbitol	1605
Glutamic acid	30125	Spectinomycin dihydrochloride	492
Glycine	30144	Threonine	30138
Guanethidine sulfate	610	Treosulfan	9126
Heptaminol hydrochloride	1428	Urea	3078
Hexandicarboxylic acid	30151	Valine	30129
Hydroxyproline	30126	Xylitol	3119
Iloprost tromethamine	30169	•	-

Pharmaceutical substances showing a well distinct UV spectrum in alkaline solution or being insoluble in methanol

Maximum in 0.1M NaOH	A 1%, 1cm	Minimum (nm)	Pharmaceutical substances	Spectrum No.
406nm	1075	253	Amphotericin B	481
372nm, <u>268nm</u>	440	323, 250	Diosmin	30102
370nm	459	312	Tenoxicam	744
347nm,248nm	393	285, 229	Isoxicam	141
340nm, <u>266nm</u>	1320	332, 280	Oxolinic acid	998
335nm, <u>272nm</u>	1130	296	Norfloxacin	997
332nm, <u>274nm</u>	1680	299, 243	Piromidic acid	959
331nm	222	280	Mesalazine	1619
325nm	68	294	Levothyroxine sodium salt	2309
324nm	75	293	Thyroxine	2305
320nm	1015	283	Tioguanine	980
319nm	355	271	Tizanidine hydrochloride	2462
304nm	809	240	Thiobutabarbital Sodium	1546
304nm	54	290	Vancomycin hydrochloride	491
294nm	130	270	Tyrosine	3007
291nm, <u>247nm</u>	286	276, 230	Amoxicillin trihydrate	433
287nm	207	240	Capreomycin disulfate	475
280nm	580	242	Diazoxide	621
279nm	700	243	Lonazolac calcium salt	725
279nm	580	241	Orazamide	30114
265nm	8.4	260	Gallamin Triethiodide	2405
262nm	232	240	Inositol Nicotinate	024
252nm	189	236	Brallobarbital calcium salt	1533
252nm	325	234	Phenobarbital	1508
251nm	274	234	Heptabarb	1529
250nm	268	232	Cyclobarbital calcium salt	1503
250nm	303	234	Allobarbital	1527
247nm	243	231	Narcobarbital	1537
247nm	298	233	Vinylbital	1523
245nm	343	230	Methylphenobarbital	1518
244nm	318	230	Amobarbital	1513
244nm	319	233	Cyclopentobarbital	1534
243nm	333	232	Aprobarbital	1514
243nm	342	230	Barbital sodium salt	1501
243nm	328	229	Hexobarbital	1505
243nm	338	230	Pentobarbital	1520
243nm	248	230	Secobarbital sodium salt	1538
242nm	335	230	Butobarbitone	1502
242nm	205	233	Propallylonal	1509
232nm	570	225	Phenformin hydrochloride	507
230nm	746	224	Buformin hydrochloride	502
228nm	1079	224	Bemegride	001
224nm	375	<215	Allantoin	1108